

Validation documentation
*APIs & excipients, liquid/semi-solid (with analysis
certificate)*

HiperScan GmbH

December 13, 2018

Contents

Introduction	7
Context of this document	7
Criteria for the inclusion of substances	8
Validation concept	9
Model creation procedure and validation runs	9
Particularities of individual substance classes	11
Significance of testing with <i>Apo-Ident</i>	12
Conclusion	13
Explanation of terminology	14
Summary	19
Validation samples	19
Validation results	19
Validation reports	21
Acetic acid 30%	21
Acetic acid 99%	25
Acetone	31
Acorus oil	41
All-rac-alpha-tocopheryl acetate	45
Angelica root oil	51
Anise oil / Star aniseed oil	55
Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)	61
Base cream DAC	67
Basis Cordes® RK	121
Benzalkonium chloride solution 50%	131
Benzyl alcohol	137
Benzyl benzoate	141
Benzyl nicotinate	145
Bergamot oil	149
Cajeput	153
Camphor oil 20%	157
Camphor ointment	161
Camphor spirit	165
Caraway oil	171
Castor oil	177
Cedar leaf oil	191
Cedar wood oil	195
Cinnamon bark oil, artificial	199
Cinnamon bark oil, Ceylon	203
Citronella oil	207
Citrus oils	211
Clotrimazole 10% Cordes® RK	219
Clove oil	223
Coal tar spirit SR	229
Collodion 4%	235
Cooling cream / Cooling cream with rose oil	239
Coriander oil	267
D-alpha tocopherol	271
Dexpanthenol	275
Dexpanthenol concentrate 50%	287
Dextrorotary lactic acid 20%	293
Diethyl ether	297
Dimethyl sulfoxide	303
Dimeticon 100 / 350 / 1000	311

Dwarf pine oil	317
Ethanol 90% pure / denatured	323
Ethanol 96% pure / denatured	333
Ethyl acetate	345
Eucalyptus oil / Rosmary oil	349
Eugenol	359
Fennel honey (conserved)	363
Fennel oil, bitter	367
Gel Cordes®	371
Geranium oil, genuine	375
Glucose syrup	379
Glycerine ointment (conserved)	383
Glycerol >99%	387
Glycerol 85%	393
Glycolic acid 70%	413
Guaiazulene 25%, water-soluble	417
Hard fat 649	423
Honeys	427
Hydrophilic base emulsion (NRF p.25)	433
Hydrophilic ointments	443
Hydroxyethyl salicylate	527
Isopropyl alcohol 100%	531
Isopropyl alcohol 70%	557
Isopropyl myristate	579
Isopropyl palmitate	585
Juniper oil	591
Lactic acid	595
Lanette® ointment (conserved)	601
Lanolin DAB	607
Lauromacrogols	619
Lavender oil / Clary sage oil	641
Liniment ointments	651
Lotio Cordes®	675
Macrogol 8 stearate type I	679
Macrogol ointment	683
Macrogol-2-laurylether	689
Macrogols 300	693
Macrogols 400	697
Mandarin oil	705
Matricaria oil, rich in (-)-alpha-bisabolol	709
Melissa oil	713
Methyl salicylate	719
Metronidazole 25% in white petrolatum	725
Mint oil / Peppermint oil	729
Myrtle oil	737
Non-ionic hydrophil ointment	741
Octyldodecanol	749
Oleic acid	759
Oregano oil	763
Paraffin classification	767
Pasta Cordes®	851
Permethrin 25% solution	857
Petroleum ether	863
Petroleum oil, white, cohobated	877
Phytomenadione	881
Pine silvestris oil	885
Pine tar	889

Polyhexanide solution 20%	893
Polyol fatty acid ester	899
Polysorbate 80 / Macrogolglycerolricinoleat	903
Potash soap DAC	909
Propylene glycol	913
Psorisol base cream Audor®	929
Rose oil (genuine / artificial)	933
Sage oil	939
Salicylic acid 50% (paraffin)	943
Shea butter, refined	953
Sodium lactate solution 50%	959
Sodium laureth sulfate solution 27%	965
Soft ointment	969
Softisan® 378	981
Sorbitan- and glycerol monooleate	987
Sorbitan trioleate	991
Spearmint oil	995
Spirit of ether	999
Spirit soap	1003
Spruce needle oil	1007
Squalane	1011
Syrups	1015
Tea tree oil	1031
Thyme oil	1035
Triglycerides	1041
Triglycerides, medium-chain	1105
Triglycerol diisostearate	1117
Verbena oil, nature-identical	1121
Viscous base solution (NRF p. 20)	1125
Water and aqueous solutions	1137
Wines	1161
Wood oils	1167
Wool alcohol cream DAB/SR	1173
Wool fat	1201
Wormwood oil	1215
Ylang-ylang oil	1219
1,3-butanediol	1223
2-ethylhexyl laurate	1227
Appendix	1233
Additional calibration samples (<i>Type A</i>)	1233
Additional validation samples (<i>Type B</i>)	1235
Additional validation samples (<i>Type C</i>)	1243
Requirements of validation	1279
Conformity of <i>Apo-Ident</i> with the <i>European pharmacopeia</i>	1280
Addition to the two-stage models	1281
References	1283
Index	1285

Introduction

The unequivocal proof of the identity of pharmaceutical raw materials based on a monograph or traditional alternative methods is work-intensive, time-intensive and economically often no longer makes sense. Near-infrared spectroscopy (NIR) provides a new means here. It enables relatively easy, fast and nonetheless safe identity testing via the preparation and evaluation of spectra.

The analysis system *Apo-Ident* was developed specially for pharmacy use. Pharmacists have the duty to check the identity of all raw materials for extemporaneous products in their pharmacies. This normally takes place based on the monographs for the respective substances in the European Pharmacopoeia. But also NIR spectroscopy is described in the European Pharmacopoeia as an identification method which, as opposed to the methods incorporated in the respective monographs, is approved for testing [1] (quotation translated)

with the prerequisite that the same results (“namely the determination of the identity” [2]), (quotation translated) are achieved as with the described methods and instruments.

The *Apo-Ident* analysis system serves to identify raw materials for prescriptions in the defined manner according to *ApBetrO* [Pharmacies Rules and Regulations] §§ 6 and 11 at pharmacies (NIR spectroscopy as an alternative testing method). *Apo-Ident* consists of three components:

- An *NIR spectrometer*, which records the spectra of non-preprocessed raw materials in a measuring glass in diffuse reflection or transfection.
- The *QuickStep* spectroscopy software controls the instrument and records the spectra and user inputs via a pharmacy-specific software plug-in. It also generates the test protocol for documentation of the testing and storage of the printout to be signed in the pharmacy.
- The software module *IdentModule* incorporates *reference databases*. The spectra from the *QuickStep* software are presented to it for evaluation.

NIR spectroscopy is a very powerful analytical method. It is also able to establish the identity of several chemical compounds and mixes in as far as an appropriate database (technically correct: a [chemometric model](#)) was created. Identity testing with *Apo-Ident* is a very safe, very fast and easy to operate analytical method for testing a large number of raw materials.

Context of this document

The suitability of the instrument, method and database is proven as follows:

- *NIR spectroscopy as a method for identity testing*: The *Ph. Eur.* [3] describes NIR spectroscopy in *Section 2.2.40* as an analytical method which is also suitable for the identification of raw materials. Therefore, validation of the method as such is not necessary.
- *Performance of the instrument*: The *Ph. Eur.* [3] furthermore describes the apparatus and the testing of its performance in *Section 2.2.40*. The document *Erfüllung von 2.2.40 Ph. Eur. durch Apo-Ident* [4] compares the implementation by *Apo-Ident* with this monograph in order to prove that *Apo-Ident* meets the specifications of the Pharmacopoeia. Each individual instrument delivered to a pharmacy is qualified in accordance with the tests described in “*Control of Instrument Performance*”. In this test, the unit consisting of analysis instrument hardware and the *QuickStep* spectroscopy software is assessed. The result is documented in a test protocol which is kept at the pharmacy.
- *Validation of the database* is documented separately for each substance class. The report at hand documents the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)*.

The *Arbeitsgemeinschaft der Pharmazieräte Deutschlands (APD)* [Working Group of German Pharmacy Inspectors] has clarified the following in its resolution dated October 16, 2013 ([5], quotation translated):

NIR is a testing method incorporated in the Pharmacopoeia. The testing quality depends on the quality of the database stored. The APD views the use of NIR instruments in case of ensured validation of the databases used in conjunction with it as one of several options for identity testing.

The APD defined more precisely ([6], quotation translated) on October 1, 2014:

The use of near-infrared is a recognised testing method according to Ph. Eur. 8. For the use of NIR instruments in pharmacies for testing the identity of raw materials, sufficient and verifiable validation of the instrument used is required. The quality of the database stored by the instrument manufacturer is decisive for quality. Batch-specific differences with the same original substances must be taken into account if present.

So NIR is basically suitable. The validity of the reference database is proven with the existing validation documentation.

Criteria for the inclusion of substances

This validation documentation describes the results of the validation of the reference database for the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)*. Validation documentation is created for each published version of the reference database for all substance classes incorporated.

The reference database is incorporated in the software module *IdentModule*. During identification testing with *Apo-Ident*, spectra which are used for evaluation purposes are presented to it by the *QuickStep* software. In the same manner, the *IdentModule* is presented all validation spectra successively during the validation runs for evaluation purposes. The *IdentModule* responds respectively (without taking the initial assumption into account) with the identified substance or rejects it as unknown. The correctness of this response is checked for each possible initial assumption and counted.

The results are summarised for each substance and reproduced in this document. The core statement of this validation report is that the following criteria must be fulfilled for each database entry, so that *Apo-Ident* offers verification of identity by means of NIR for the relevant substance/substance group:

- The database is exclusively generated from spectra which have been recorded by *HiperScan GmbH* on traceable samples in pharmaceutical quality.
 - The samples are procured via typical pharmacy sources (*DAC III.2.: Bezugsquellennachweis für Rezepturbestandteile [reference source for prescription components]* [7]).
 - A valid manufacturer's certificate exists (content, purity and identity of the batch).
 - The identity was confirmed by a certified test laboratory or *HiperScan GmbH*.
- Each version of the reference database (every update) is validated in-full.
 - Calibration spectra (*Type A*), other spectra recorded under the control of *HiperScan GmbH* (*Type B*), and spectra from the field (*Type C*) are presented to the *IdentModule* for evaluation in three separately evaluated validation runs.
 - Here, no single *false positive* result may arise.
 - Here, the various substance classes are also tested for reciprocal rejection, where this is objectively justified (see *Summary* section).
- In the validation with spectra recorded under the control of *HiperScan GmbH*, spectra of at least one independent sample must be considered, i.e. spectra from a batch whose spectra have not been used for the generation of the database. In addition, the set of *Type A* and *Type B* spectra must originate from at least three different batches.
- Spectra of additional substances may be used for the generation of the database even though they will not be offered for identification by this database. The purpose is the reliable distinction from these substances.

- Any positive result of *Apo-Ident* confirms the identity of the substance/substance group and distinguishes it from all other substances of the database. In the case of substance groups, the result is ambiguous: Distinction from all substances not belonging to the group is proven. The substance is identified as a member of this group. However, within the substance group, it is not possible to reliably classify which substance has been tested.
- The criteria for clear identifiability are a [specificity](#) of 100 % ([true negative rate](#)) and a minimum distance in the distance matrix. See 2. d) under [Model creation procedure and validation runs](#).

Validation concept

[Chemometrics](#) is a statistical technique for the extraction of relevant chemical information from spectra. In mathematics, this method is described as *multivariate data analysis*. Chemometrics proceeds here as follows:

1. Collection of spectra for the *calibration sample*. The results (identities) of the calibration sample must be known. The calibration samples must be representative for the samples which are to be evaluated later. Therefore, they must take the various possible (physical) compositions into account. (Therefore, sourcing calibration samples for NIR from the specialist trade is superior to the use of CRS reference substances.)
2. The first mathematical step is *calibration*. Here, the [chemometric model](#) is calculated from the *calibration sample spectra* ([reference spectra](#)) and limits as well as some parameters are stipulated. The chemometric model is used later to calculate the analysis result (*prediction*).
3. Collection of further spectra for the *validation sample* which should be independent of the *calibration sample*. The results (identities) of the *validation sample* must also be known. The textbook suggests a random sample with a normal scope of 25 % to 50 % of the *calibration sample* [8].
4. The second technical data step is *validation*. Here, the [chemometric model](#) created is evaluated based on the spectra of the *validation samples*. As validation parameters for the identification, the *Ph. Eur. Section 2.2.40* [3] specifies the [specificity](#) and [robustness](#).

The validation step according to the textbook has the target of estimating the performance capability of the model created based on a random sample. In order to achieve the best possible precision, attention is paid to the calibration sample. In the field of pharmaceuticals, the safety of the method has priority. In order to be able to [validate](#) the model within the regulatory scope, the validation step must include probative force. For this purpose, the validation sample must be *representative and complete* in order to enable the testing of all cases.

A *sufficient number of batches* must be secured for validation because validation finally proves whether the number of batches in calibration suffices.

Each substance is validated individually. The validation results are documented per substance in this document. Moreover, the documents show how many and which batches have been used for creating the model or model validation.

At least one certificate is taken in for each substance from an accredited test laboratory for the independent testing of identity of the sample. The identification number of the corresponding test certificate is listed in the report, enabling traceability of a substance tested according to the monographs in the Pharmacopoeia.

Model creation procedure and validation runs

The safety of the [chemometric models](#) is guaranteed by several measures during model creation, of which the validation step is the final one. Normally, the procedure is as follows. It is in particular valid for the active pharmaceutical ingredients (APIs) *solid API excipients, liquid/semi-solid API excipients (with a test certificate), narcotics - solid medicinal substances and drugs*. If, for individual substance classes, variations are required, they are depicted in the section [Particularities of individual substance classes](#).

1. Collecting the reference spectra (calibration sample)

- a) Procurement of the samples from the same sources from which pharmacies source their raw material for compounding (Caelo, Fagron, Euro-OTC, . . . , see also *DAC III.2. Bezugsquellen-nachweis für Rezepturbestandteile* [Sources of supply for compounding] [7]).
- b) Testing the suitability according to *ApBetrO* [Pharmacies Rules and Regulations] §§ 6, 11, that is to say the availability of a valid manufacturer certificates via identity, purity and contents of the batch.
- c) Recording standard 40 spectra of the sample in different positions, as a standard on four instruments. Here, handling and presentation of the samples as later in the pharmacy.
- d) Visual checking for anomalies in the spectra. In case of indications of measurement errors, measurement must be repeated. If a signature is missing in the spectrum, the substance may be excluded from the start as not promising (the spectra are nonetheless entered in the database validation as independent *Type B* spectra).
- e) Testing identity. For each substance, a certificate of correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the respective following substance page of this validation documentation the *Mahalanobis distance* to this reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance. Such samples underline the statistical spread of the original reference substance, but cannot add any new characteristics of the substance.

HiperScan GmbH cooperates with some suppliers as follows: the raw materials supplier takes a sufficiently large sample in his incoming goods area so that part of it can be used for recording the NIR spectra. The remainder of the sample goes to analytics for market approval. The manufacturer's batch certificate emerges from these identity, contents and purity tests, which consequently also prove the correct identity of the NIR reference sample. Therefore, the NIR spectra are suitable for structuring the database (*Type A*) and can also optionally be used for validation (*Type B*). The samples which this applies to are marked in the validation report with a footnote.

- f) If the identity of the new sample is proven, it is declared as a reference sample and the spectra are approved for structuring the database.

2. Generating the chemometric models (calibration)

- a) Determination of the transformation matrix from the reference spectra using variance maximisation [8, 9]. (All reference spectra are always included, even if only a few spectra are added for an update.) All reference spectra receive the same data pre-treatment, which is also later applied to all measurement spectra in the field (in the pharmacies).
- b) Checking that the number of principal components used is still sufficient.
- c) Calculating the limits for each substance from the spread of the reference spectra. The calculation regulation is identical for each substance in a substance class.
- d) Checking the distances between the limits of the separable substances: the distance matrix incorporates the *Mahalanobis distances* from each substance to all the others. If a distance is less than 10, both two substances must be declared as inseparable or *both two* must be removed from the database. (The spectra of removed substances remain as independent *Type B* spectra in the validation.)
- e) Testing the model based on the reference spectra. If *false positive* results arise, you must proceed as in the case of *Mahalanobis distances* which are too small.
- f) If chemometric models are available for all substance classes which meet both criteria (distance matrix and no *false positives*), they are joined together with the evaluation algorithms as an *IdentModule* and encrypted. This unit can no longer be changed. Its overall function is tested by validation.

3. Set of validation spectra (validation samples)

The following is provided for validation:

- a) *Type A*: The reference spectra = calibration spectra from which the database was generated. These also include spectra from substances which the *chemometric model* should not identify, but were also recorded during generation in order to increase selectivity. (As a result, the model “learns” to differentiate from other substances which are actually unknown to it.)
- b) *Type B*: Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database. These also include reference spectra of other substance classes, and spectra that are not used as reference spectra. Samples are considered to be independent, if they originate from a batch, of which no spectra have been used for the generation of the database. (Up to *IdentModule 2018-01*, samples were still considered independent if the sampling was done independently, i.e. if they originate from another sales container.)
- c) *Type C*: Spectra from the field, which have not been recorded under the control of *HiperScan GmbH* and have not been used for the generation of the database. The spectra include not only substances of the substance class to be tested, but also substances from other classes.

All manufacturers’ batches from which spectra are used for the validation are listed by substance in this document: for substances included in the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)* in the respective validation reports; otherwise in attachments *A*, *B* and *C*.

Furthermore remains valid: validation spectra may only be removed if a error in the spectrum can be proven. Here, the spectra are not deleted, but instead placed on a *blacklist* incorporating the reason, date and initials in the commentary.

The section *Particularities of individual substance classes* treats the other substance classes from which *Type B* and *Type C* spectra are cited for validation purposes.

4. Validation runs and approval

- a) Validation spectra are transferred holistically to the *IdentModule* for evaluation in the same way as the spectroscopy software *QuickStep* transfers measured spectra.
- b) Following the provision of each spectrum, the *IdentModule* responds as to whether it has recognised a substance and which substance was recognised.
- c) The correctness of this response is checked for each possible initial assumption (each measurable substance with the substance class) and counted according to *true negative*, *false negative*, *true positive* and *false positive*. These figures are provided for each substance and additionally in the section *Summary*, separated according to types *A*, *B* and *C*.
- d) No *false positive* results whatsoever are permissible.
- e) If the criterion is also met for all substance classes, the *IdentModule* is approved.

Particularities of individual substance classes

Basically, *HiperScan GmbH* procures and tests the manufacturer’s certificate for the batch, commissions external testing of the identity of the sample or carries it out independently and stores the certificates. As described, this process is established for the Pharmacopoeia substances, that is to say for substance classes **APIs & excipients, solid**, **APIs & excipients, liquid/semi-solid (with analysis certificate)**, **Narcotic substances, liquid/semi-solid** and **Drugs**. Therefore, *HiperScan GmbH* is able to furnish proof of the identity of the reference samples. In case of manufacturer-specific substance classes and others, individual steps are organised differently in-part:

The substance class **APIs & excipients, liquid/semi-solid (other)** (often described as cosmetics) incorporates substances for which no specification of the requirements of the pharmaceutical quality is determined, neither in a Pharmacopoeia monograph, a DAC/NRF monograph nor via a manufacturer’s specification. Consequently, neither the identity nor contents can be tested independently. No certificates whatsoever exist for the reference samples. So here, merely the matching of the sample with former samples of this product is established and confusion with the other substances

is ruled-out. (If the manufacturer of such a substance prepares a specification, determines testing methods and provides manufacturer's certificates in accordance with *ApBetrO* [Pharmacies Rules and Regulations] §§ 6, 11, *HiperScan GmbH* can assign the substance to the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)* again in the future).

Substance class **HCK – nutritional supplements (Hepart)** contains the HCK micro-nutrients from the Swiss company *Hepart AG*. *HiperScan GmbH* receives the reference samples directly from the manufacturer. For each reference sample, *HiperScan GmbH* also receives manufacturer's certificates and keeps these. New checking of the identity of the reference sample is not carried out by *HiperScan GmbH*. The identity of the reference samples is therefore documented by *Hepart AG*. The spectra of all batches provided by *Hepart AG* are recorded by *HiperScan GmbH* and entered in the database.

All the manufacturer's batches are used for the generation and validation of the substance class *HCK*. The expected variation is also represented in the generation and validation if there are less than three batches.

Also, for the substance class **PhytoComm** (TCM-Granulated herbal extracts of the manufacturer *PhytoComm*) spectra for all useable batches are recorded by *HiperScan GmbH* and entered in the database. The supplier organises the respective tests themselves and keeps the test certificates.

A new evaluation option was created for the class *PhytoComm* with the update 2016-01. As the risks are considerably fewer than those from chemical agents, the pharmacist can specify a reasonable criterion for the *specificity* in accordance with internal risk estimation. The database for this is created without taking safety distances into account and no criterion is determined in advance for the *specificity*. Instead, the *specificity* for testing the identity with this concrete substance is calculated in the validation for each substance and provided with the measurement result. The pharmacist then judges himself whether this safety is reasonable with regard to the risk of the substance.

Additionally, a statistical forecast is provided for the *specificity* which is determined according to the *Rule of Three* [10, 11]. For this forecast, it is assumed that there would have been three wrong results more and is provided with a lower limit for *specificity*. This value has a special meaning if a *specificity* of 100 % is achieved for a substance during validation. In this case, the lower limit allows conclusions regarding the scale of existing safety for which with an endless number of validation spectra a value of less than 100 % is to be assumed.

If, for example 14000 spectra not belonging to the substance are presented and no *false positive* classification is made, a hypothetical number of three *false positive* results is assumed (*Rule of Three* [10, 11]) and the *specificity* is defined with 100.0000 % (> 99.9786 %). Here, it applies that the higher the number of validation spectra which form the statistical basis, the better the *specificity* calculated via the lower *specificity* limit will be approximated.

The positive result of the identity test using *Apo-Ident* establishes that the sample spectrum is in accordance with a batch of the specified granulate from the supplier *PhytoComm*, whereby all useable batches from the supplier are known.

The *PhytoComm* class can only confirm the identity of batches that have been used for the generation of the database. As a consequence there cannot be any validation spectra of other batches. Therefore, the criterion reads that two samples (from different sales containers) from each batch must exist, one for the structure of the database (*Type A*) and one for the validation (*Type B*).

Significance of testing with *Apo-Ident*

The analysis result is determined using sophisticated statistical methods according to state-of-the-art science and technology. Chemical and pharmaceutical knowledge is applied for the selection of the samples from which the calibration spectra and validation spectra are recorded. Otherwise it does not influence the further steps of model creation.

Verbally, the statement of the analysis result can be expressed as follows. Here "*the spectra match*" means that the criteria *Mahalanobis distance*, *outlier analysis* and *correlation* are met as shown in *Erfüllung von 2.2.40 Ph. Eur. durch Apo-Ident* [4]. "The spectra do not match", on the other hand, means that at least the criterion *Mahalanobis distance* is not met.

The positive analysis result "*was identified as ...*" is very meaningful because both the quantity of substances to be taken into account and the number of underlying samples is very comprehensive.

1. The spectrum of the sample measured matches spectra of the defined substance.

2. The spectrum of the sample measured does not match any spectrum of any other substance in this substance class. Therefore, all other substances can be clearly ruled-out.
3. As the spectra from other substance classes were used for validation, it is proven that no spectrum of one of these other substances matches the defined substance. (All substance classes with which a spectrum comparison is possible and makes sense are used for validation. This is documented for each substance class in the section [Summary](#).)
4. If the defined substance belongs to a [substance group](#) which in itself is not clearly separable with *Apo-Ident*, matching with the spectra of one or several substances in this group is confirmed. Which of these substances it actually is cannot be determined clearly. All other substances are excluded analogous to [2](#) and [3](#).

On the other hand, a negative analysis result “*was not identified as ...*” means:

1. The substance offered could not be recognised based on the spectrum of this sample.
2. The identity of this sample is not confirmed.
3. Testing must be repeated in accordance with the specifications of the Pharmacopoeia.

Conclusion

NIR spectroscopy is a testing method incorporated in the Pharmacopoeia. In case of successful database validation, it is a possible method for identity testing [\[5\]](#). *Apo-Ident* meets the criteria of the *European Pharmacopoeia* as a near-infrared spectrometer and proves the validity of the reference database with the existing validation documentation. This means that *Apo-Ident* can be used as an alternative testing method for testing raw materials at pharmacies.

Explanation of terminology

The following section serves to explain or define specialist terminology which is required in order to understand this document. If necessary, definitions for the analysis system *Apo-Ident* are defined more precisely.

The term database is used in this document exactly as in the *Ph. Eur. Section 2.2.40* [3] synonymous with [chemometric model](#). In order to differentiate the databases which are relatively independent of each other, *HiperScan GmbH* frequently also uses the term [substance class](#) (primarily in the plural). On the other hand, the spectra used to structure the database are termed spectrum collection and not database.

Substance classes are units of the organisational structure of the *IdentModule*. The substance classes are substance [databases](#) which are also broadly independently subscribable. On the one hand, the liquid and semi-solid substances are separated from the solid powders because they are measured against different references and therefore the spectra cannot be compared. On the other hand, for example the Pharmacopoeia substances are kept separated from the manufacturer-specific database *PhytoComm* for TCM (traditional Chinese medicine) raw materials.

The individual substance classes need only be limited against each other in-part. Often, no risk of confusion exists because they can only be procured from different sources. On the other hand, in several cases we handle substances which need not be distinguished. For example, em Huang Qi granulate from the company *PhytoComm* neither needs to be delimited from *Huang Qi* granulate from the company *HerbaSinica* nor is matching required. Respectively one single [chemometric model](#) is behind a substance class. (Even if several reciprocally secured chemometric models would be permissible.) The terms *substance class*, *chemometric model* and *databases* are mostly used here as synonyms.

A substance group respectively summarises all the substances within a [substance class](#) which cannot safely be distinguished from one another based on their NIR spectra. However, all the other substances in the database can be excluded.

The formation of subgroups is mentioned in the *Ph. Eur. Section 2.2.40* [3]. In this manner, technical restrictions in case of extensive databases can be avoided and it is possible to prepare individual subgroups with different spectrum pre-treatment. Validation of the subgroups against each other is required. *HiperScan GmbH* has solved these technical restrictions and doesn't use any subgroups within a substance class any longer.

Principal component analysis (PCA) [8, 9] is a multivariate statistics process or multivariate data analysis. It serves to structure, simplify and illustrate comprehensive data records by describing a large number of statistical variables by describing a lower number of linear combinations (the *principal components*) which are as significant as possible. In the *Apo-Ident IdentModule*, *PCA* is used to evaluate the recorded spectrum data (corresponding with *Ph. Eur. 2.2.40* [3]).

The term validation is defined in both relevant contexts here with different (even if related) meanings.

Within the sense of the expert discipline of *chemometrics*, validation is a process step when creating a [chemometric model](#): after a transformation matrix, limits and various parameters have been calculated or determined from a set of reference spectra during the course of the calibration step [8, 9], the validation step determines the performance capability of the model (selectivity, precision, ...) based on the validation spectra. Normally, random sampling is planned here. In order for the validation to gain strength of proof, the validation spectrum set must be selected with an appropriately wide scope (*representative* and *complete*). The terms *validation run* and *validation step* always actually mean the process step in this sense.

In the regulatory sense (of pharmaceutical production), validation is the documented proof that a process or system meets the previously specified requirements reproducibly when applied practi-

cally. In this sense, the *Apo-Ident* databases only become validated databases with the validation documentation, which this document is part of.

The *European Pharmacopeia* uses the term validation in *Section 2.2.40* within the sense of the specialist discipline of *chemometrics* [3].

The robustness of a process is the property of only being influenced by environmental fluctuations (e.g. temperature or humidity) a little. A method is robust if the environmental conditions do not or hardly falsify the final result.

The specificity of a classification (of a [chemometric model](#)) is the [true negative rate](#).

The recognition rate (also sensitivity) is the [true positive rate](#). It defines in how many percent of cases a correctly set up substance is actually confirmed.

The true negative rate describes the share of spectra correctly classified as non-identity during validation. This is equivalent to correct classification. It means that a substance *A* within identity checking as substance *B* is judged as “*not identified*”. The *true negative rate* is equivalent to the conditional frequency

$$h(\text{rejected}|\text{genuinely no identity}) = \frac{r_n}{r_n + f_p}$$

with r_n as the total number of *true negative* classifications and f_p as the total number of *false positive* classifications. For successful validation of an *IdentModule*, all spectra presented belonging to this category must be classified as *not in accordance*.

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of this number. The weight of each spectrum of a substance/substance group *i* therefore results as

$$w_i = \frac{1}{n_i}$$

with n_i number of spectra in this substance/substance group. This weighting ensures that the overall result cannot be enhanced by adding especially large numbers of spectra from easily separable substances.

The true positive rate describes the share of spectra correctly classified as identity during validation. This is equivalent to correct classification. It means that a substance *A* within identity checking as substance *A* is judged as “*identified*”. The *true positive rate* is equivalent to the conditional frequency

$$h(\text{identified}|\text{genuine identity}) = \frac{r_p}{r_p + f_n}$$

with r_p as the total number of *true positive* classifications and f_n as the total number of *false negative* classifications. The *true positive rate* is a measure for the recognition rate of the validated *Apo-Ident IdentModule*.

In order to ensure that each substance is received with the same weight, the spectra are weighted as described for the [true negative rate](#).

The true negative result describes a spectrum correctly classified as non-identity during validation. It is equivalent to correct classification. It means that a substance *A* within identity checking as substance *B* is judged as “*not identified*”.

The false positive result describes a spectrum falsely classified as non-identity during validation. This is the most critical type of possible false classification. It means that a substance *A* within identity checking as substance *B* is judged as “*identified*”. For successful validation of an *IdentModule*, a number of false positive events of zero are demanded for all spectra entering the validation. The exception to this restriction is the class of TCM granulates from the company *PhytoComm* as described under [Particularities of individual substance classes](#).

The true positive result describes a spectrum correctly classified as identity during validation. It is equivalent to correct classification. It means that a substance *A* within identity checking as substance *A* is judged as “*identified*”.

The false negative result describes a spectrum falsely classified as non-identity during validation. It is equivalent to false classification. It means that a substance *A* within identity checking as substance *A* is judged as “*not identified*”.

The ‘Rule of Three’ says that with a probability of 95 % the next random sample of the same size no more than three false results are to be expected if no false result existed in the existing random sample [10, 11].

The *specificity* and *recognition rate* are determined both globally and from the validation runs for all substances. The information is supplemented with the hypothetical value if there had been three false results more. The percent information is provided in parentheses with the “greater than” symbol ‘>’, e.g. *specificity* 100.000 % (>99.983 %) if 17 567 false spectra have been presented without one single *false positive* result. The larger the statistical basis, the lower the influence of the hypothetical false results.

The Mahalanobis distance is a distance measure between two points in *n*-dimensional vector space. Here, the respective direction component of the distance to *standard deviation* [12] of an *n*-dimensional distribution is standardised. In case of the *principal component analysis* [8, 9] this standardisation relates to the distribution of the respective calibration data set for a classification (substance/substance group) in the *principal component space* [8]. The *Mahalanobis distance* of a point (mapping of a spectrum) $\vec{y}$ in the *n*-dimensional principal component space to the expected value of an *n*-dimensional distribution $\mathbf{X}$ then results as

$$d(\mathbf{X}, \vec{y}) = \sqrt{(\bar{\mathbf{X}} - \vec{y})^T \mathbf{S}^{-1} (\bar{\mathbf{X}} - \vec{y})} \quad \text{with} \quad \mathbf{X} \in \mathbb{R}^{m \times n}, \vec{y} \in \mathbb{R}^m$$

[13]. Here, *m* is equivalent to the number of principal components used (dimension of the principal component space) and *n* the number of measurements existing in the calibration data set (spectra). $\bar{\mathbf{X}}$ is the expected value of the resulting distribution for the calibration data set (the average value of *n* measurements received). $\mathbf{S}^{-1}$ is the inverse covariance matrix [12] for distribution $\mathbf{X}$.

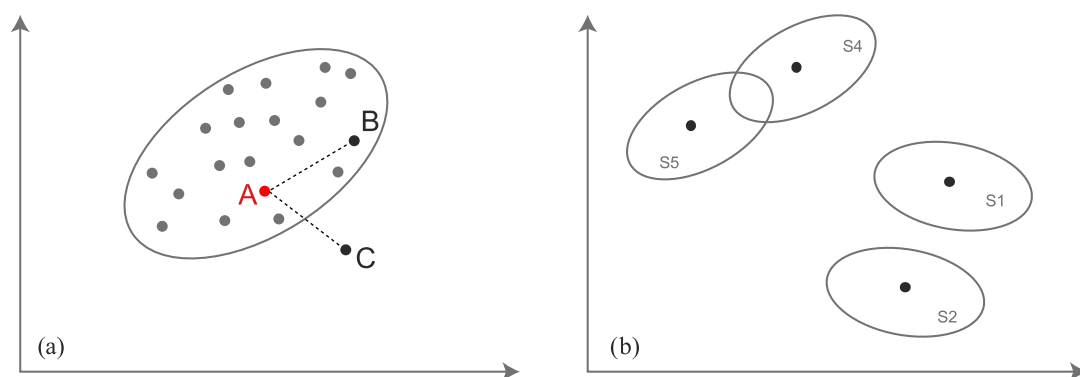


Figure 1: (a) The *Mahalanobis distance* from *A* to *B* is less than from *A* to *C*. However, the Euclidean distances are identical. (b) The *Mahalanobis distance* between the two measurement series *S4* and *S5* is smaller than between *S1* and *S2*. However, the Euclidean distances are identical.

The *Mahalanobis distance* offers advantages compared to the *Euclidean distance*: For the calculation of the distance it takes the statistical properties of a data point distribution (measurement series), i.e. average value, variance and covariance of the data points [14] into account. The *Mahalanobis distance* is applied while creating the reference database for evaluating the spectra from different samples of a substance.

A chemometric model is a classifier based on statistical methods [8, 9]. Through the respective algorithm used (e.g. *Principal Component Analysis*, *Cluster Analysis*), a maximum of chemical information is extracted from measurement data. Here, systematical or physical disturbances are eliminated using appropriate data pre-processing [15, 16].

At several places in this document, in order to simplify understanding, the term **database** is used instead of *chemometric model* – in the same manner as in the *Ph. Eur. Section 2.2.40* [3].

A sample (with its own sample ID) refers to substance in a sales container. Repeated sampling from the same sales container is listed under the same sample ID. (The suffix “SI” is not part of the sample ID.) Several samples may originate from the same batch. Samples are called “independent”, if they originate from a batch, of which no spectra have been used for the generation of the database. (Up to *IdentModule 2018-01*, samples from different sales packages were considered to be independent.) The information above the list of validation spectra now includes also the number of batches that deliver independent samples for the validation (for both *Type B* and *Type C*).

In case a supplier takes a sample for testing from its incoming goods and splits it to multiple laboratory containers, the substance in all laboratory containers will still be ascribed to the same sample. *HiperScan GmbH* only uses one of the subsamples.

Reference samples are used to structure the database. The *reference spectra* originate from these samples. In chemometric technical jargon you would normally say: For *calibration*, a **chemometric model** is generated from the *calibration spectra* recorded from the *calibration samples*, whose quality is subsequently assessed in **validation**.

Reference samples are procured via typical pharmacy sources. Their identity is tested. The *reference spectra* are recorded by *HiperScan GmbH*. The documentation also includes the manufacturer’s name and batch number.

Reference samples are clearly identified by a sample ID. Samples without sample ID may not be used as *reference samples*.

Summary

A total of 220 770 spectra from 29 736 different batches for a total of 441 substances were used to validate the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)*.

Validation samples

The validation samples can be categorised as follows:

Type A Calibration spectra. These are the spectra used to generate the chemometric model. They were recorded by *HiperScan GmbH*. Detailed information regarding the batches or samples can be found in the following validation reports under *calibration samples* and under *Type A*. Further information is listed in [Appendix A](#).

Substance class	Substances	Batches	Spectra
APIs & excipients, liquid/semi-solid (with analysis certificate)	212	840	35 544

From category *A* a total of 35 544 spectra from 840 batches for a total of 212 substances were taken into account for validation.

Type B Spectra from independent samples which are not included in database generation. These spectra were recorded by *HiperScan GmbH*. Detailed information regarding the batches or samples can be found in the following validation reports in the section *Type B* or in [Appendix B](#).

Substance class	Substances	Batches	Spectra
APIs & excipients, liquid/semi-solid (with analysis certificate)	437	2191	93 874

From category *B* a total of 93 874 spectra from 2191 batches for a total of 437 substances were taken into account for validation.

Type C Spectra from independent samples which are not included in database generation. *Apo-Ident* customers carried out the measurements. Detailed information regarding the batches or samples can be found in the following validation reports in the section *Type C* or in [Appendix C](#).

Substance class	Substances	Batches	Spectra
APIs & excipients, liquid/semi-solid (with analysis certificate)	334	27 739	91 352

From category *C* a total of 91 352 spectra from 27 739 batches for a total of 334 substances were taken into account for validation.

Validation results

The validation runs checked whether all substances/substance groups in the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, the matching of all relevant spectra of various substances with the substances/substance groups in the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)* was checked and the correctness of the results was evaluated. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	35 544	0	4 762 896
Type B	0	76 932	1055	12 591 589
Type C	0	83 374	2570	12 245 056

All substances/substance groups in the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)* can be clearly distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.000 00 % (> 99.999 78 %)	100.000 00 % (> 99.970 45 %)
Type B	100.000 00 % (> 99.999 82 %)	97.157 12 % (> 97.144 15 %)
Type C	100.000 00 % (> 99.997 48 %)	94.434 79 % (> 94.282 94 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Acetic acid 30%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30487-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Acetic acid 30%; Acidum aceticum 30 %; Acidum aceticum dilut. 30%

Special notes

When selecting the *Acetic acid 30%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Acetic acid 30%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Acetic acid 30%	12376101	31131	40	20130107
Caelo	Acetic acid 30%	12376101	31318	40	20130107
Caelo	Acetic acid 30%	16079301	33426	40	1801528

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Acetic acid 30%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 424 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Acetic acid 30%*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Acetic acid 30%	82207428	30258	40
Caelo	Acetic acid 30%	11271701	30487	40
Caelo	Acetic acid 30%	12376101	31318 [†]	20
Caelo	Acetic acid 30%	172744	33889	40
Caelo	Acetic acid 30%	18123201	34236	40

- 93 694 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 66 spectra from 53 *Apo-Ident* customers from 10 batches from the substance/substance group *Acetic acid 30%*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Acetic acid 30%	14208501	24
Caelo	Acetic acid 30%	16079301	10
Caelo	Acetic acid 30%	12376101	14
GEHE am 21.05.2014	Acetic acid 30%	12376101	1
Caesar & Loretz GmbH	Acetic acid 30%	14208501	1
Caelo	Acetic acid 30%	1428501	1
Caelo	Acetic acid 30%	14058001	1
Bombastus	Acetic acid 30%	16079301	1
Cealo/Noweda	Acetic acid 30%	16079301	2
Caesar&Loretz	Acetic acid 30%	16079301	1
Kehr	Acetic acid 30%	11271701	1
Caelo	Acetic acid 30%	11271701	3
Phoenix	Acetic acid 30%	11271702	1
Caelo	Acetic acid 30%	11271702	2
Caelo	Acetic acid 30%	10313102	1
Caelo	Acetic acid 30%	12376102	1
Phönix	Acetic acid 30%		1

- 91 286 spectra from 501 *Apo-Ident* customers from a total of 27 000 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Acetic acid 30%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Acetic acid 30%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	165	15	93 694
Type C	0	61	5	91 286

The substance/substance group *Acetic acid 30%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	91.6667 % (> 90.0000 %)
Type C	100.0000 % (> 99.6587 %)	92.4242 % (> 87.8788 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31131	33426	13.63	47.91
31318	33426	4.54	47.03
33426	33426	0.00	43.36

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Acetic acid 99%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30486-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Acetic acid 99%; Acidum aceticum glaciale

Special notes

When selecting the *Acetic acid 99%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Acetic acid 99%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Acetic acid 99%	13013101	31132	39	20130114
Caelo	Acetic acid 99%	13013101	31319	40	20130114
Caelo	Acetic acid 99%	16238401	33423	40	1801527

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 119 spectra of 3 reference samples from the substance/substance group *Acetic acid 99%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 425 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 4 reference samples of the substance/substance group *Acetic acid 99%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Acetic acid 99%	11186701	30486	80
Caelo	Acetic acid 99%	13013101	31319 [†]	20
Hedinger	Acetic acid 99%	223083B	33428	40
Hedinger	Acetic acid 99%	223083B	33429	40

- 93 694 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 105 spectra from 81 *Apo-Ident* customers from 31 batches from the substance/substance group *Acetic acid 99%*.
- These include spectra of independent samples from 29 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Acetic acid 99%	13218301	18
Devisa/Noweda	Acetic acid 99%	13218301	1
Caelo	Acetic acid 99%	13218501	2
Fishar/AHD	Acetic acid 99%	7004111	1
Fischar	Acetic acid 99%	7004051	1
Caelo	Acetic acid 99%	16238401	7
Caelo	Acetic acid 99%	11186701	1
Caelo	Acetic acid 99%	10135706	1
Caelo	Acetic acid 99%	11349901	5
Sanacorp	Acetic acid 99%	11349901	1
Caelo	Acetic acid 99%	13218401	1
Hedinger	Acetic acid 99%	223018B	1
Caelo	Acetic acid 99%	15280701	9
Caelo	Acetic acid 99%	13013101	6
Sanacorp	Acetic acid 99%	13218301	1
Hedinger/NOWEDA	Acetic acid 99%	223012B	1
Hedinger/Jenne	Acetic acid 99%	223044B	1
Hedinger	Acetic acid 99%	223052B	1
Hedinger	Acetic acid 99%	5113M-01010	1
Fischar/Gehe	Acetic acid 99%	7003053	1
Phoenix	Acetic acid 99%	7003053	1
Otto Fischar / Phoenix	Acetic acid 99%	7003053	1
Otto Fischar GmbH	Acetic acid 99%	7003053	1
Fischar	Acetic acid 99%	7003053	1
otto-fischar/gehe	Acetic acid 99%	7004051	2
Fischer/Phönix	Acetic acid 99%	7004111	1
Fischar/Phönix	Acetic acid 99%	7004111	1
otto fischar	Acetic acid 99%	7004111	1
Fischar	Acetic acid 99%	7005091	2
ACM	Acetic acid 99%	KNC1163J14	1
Caelo	Acetic acid 99%	15280702	6
Caelo	Acetic acid 99%	17200401	1
fischar	Acetic acid 99%	180406E/7006091	1
Hedinger	Acetic acid 99%	223055B	1
Hedinger	Acetic acid 99%	223083B	1
Hedinger / Kehr	Acetic acid 99%	223079B	1
Fischar Otto / Sanacorp	Acetic acid 99%	7005091	1
Fischar/Kehr	Acetic acid 99%	7005091	1
Fischar/Sanacorp	Acetic acid 99%	7006091	2
Fischar/Noweda	Acetic acid 99%	7006091	1
GEHE	Acetic acid 99%	7006091	2
Bombastus	Acetic acid 99%	7006091	1
Fischar	Acetic acid 99%	90212/1	1
Hedinger/AHD	Acetic acid 99%	2230128	1
Fischar/Phoenix	Acetic acid 99%	7002031	1
Fischar/Gehe	Acetic acid 99%	7002031	1
Otto Fischer GmbH	Acetic acid 99%	7002031	1
Anzag	Acetic acid 99%	7002101	1
ph.eur /phönix	Acetic acid 99%	7002101	1
Otto Fischar/ Phoenix	Acetic acid 99%	7002101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Acetic acid 99%	7002101	1
Fischar/Phoenix	Acetic acid 99%	7002101	1
Caelo	Acetic acid 99%	7002101	1
Fischar / Anzag	Acetic acid 99%	7002031	1
Hedinger/Gehe	Acetic acid 99%	223007B	1

- 91 247 spectra from 501 *Apo-Ident* customers from a total of 26 978 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Acetic acid 99%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Acetic acid 99%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	119	0	35 425
Type B	0	180	0	93 694
Type C	0	105	0	91 247

The substance/substance group *Acetic acid 99%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 94.9580 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.6667 %)
Type C	100.0000 % (> 99.6586 %)	100.0000 % (> 94.2857 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31132	33423	10.78	92.67
31319	33423	6.54	85.75
33423	33423	0.00	87.24

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Acetone**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30539-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Acetone; Acetinum

Special notes

When selecting the *Acetone* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Acetone*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Hedinger	Acetone	011215B	31575	40	1505477
Otto Fischar	Acetone	4203051	31580	40	20130513
Hedinger	Acetone	011242B	31793	40	20150703
Otto Fischar	Acetone	4205031	31796	40	20150305

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Acetone*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 239 spectra from 8 reference samples of the substance/substance group *Acetone*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Hedinger	Acetone	011152B	30539	40
Otto Fischar	Acetone	zu erfassen	30564	40
Hedinger	Acetone	011215B	31575 [†]	19
Otto Fischar	Acetone	4203051	31580 [†]	20
Hedinger	Acetone	011242B	31793 [†]	20
Otto Fischar	Acetone	4205031	31796 [†]	20
Hedinger	Acetone	011236B	32920	40
Hedinger	Acetone	011297B	33946	40

- 93 635 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 471 spectra from 178 *Apo-Ident* customers from 181 batches from the substance/substance group *Acetone*.
- These include spectra of independent samples from 175 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Hedinger	Acetone	11065B	1
Hedinger	Acetone	11134B	1
Hedinger	Acetone	11139b	1
Hedinger	Acetone	11139B	1
Hedinger	Acetone	11143B	3
Anzag/Dr. Seeger	Acetone	111488	1
Caelo	Acetone	11148B	2
Hedinger, Heiligenwiese...	Acetone	11148B	1
Hedinger	Acetone	11148B	1
Hedinger/Noweda	Acetone	111528	1
Hedinger	Acetone	11152B	1
Caelo	Acetone	11152b	1
Hedinger	Acetone	111588	1
Hedinger Jenne	Acetone	11158B	1
Hedinger	Acetone	11158B	4
Großhandel	Acetone	11158B	1
Phönix/ Hedinger	Acetone	11158B	1
Anzag/ Hedinger	Acetone	11158B	1
Hedinger/Gehe	Acetone	11158B	1
Hedinger/AHD	Acetone	11158B	2
Hedinger, Noweda	Acetone	11163B	1
Hedinger	Acetone	11163B	3
Hedinger	Acetone	111688	1
Hedinger , Noweda	Acetone	11168B	1
Hedinger/Phönix	Acetone	11168B	1
Großhandel	Acetone	11168B	1
Hedinger	Acetone	11168B	1
Hedinger/Noweda	Acetone	11168B	1
Hedinger/Otto	Acetone	111738	1
Caelo	Acetone	111738	1
Hedinger	Acetone	111738	1
Hedinger / Gehe	Acetone	11173B	1
Hedinger	Acetone	11173B	11
Hedinger/Anzag	Acetone	11173B	1
Anzag	Acetone	11173B	1
Hedinger Phönix	Acetone	11173B	1
Apoth.E.Bartel / Noweda	Acetone	11173B	1
Hedinger/ Noweda	Acetone	11177B	1
Hedinger	Acetone	11177B	1
Hedinger/Gehe	Acetone	11177B	1
Barthel	Acetone	11177B	1
Gehe 06.05.2013	Acetone	111828	1
Hedinger	Acetone	111828	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Acetone	11182B	2
Hedinger/Anzag	Acetone	11182B	1
Hedinger , Stuttgart	Acetone	11182B	1
Hedinger/Noweda	Acetone	11185B	1
Hedinger/Sanacorp	Acetone	11185B	2
Hedinger	Acetone	11185B	2
Caelo	Acetone	11185B	2
Hedinger/Anzag 03.06.2013	Acetone	11185B	1
Anzag/Seeger	Acetone	11185B	1
Hedinger/Phoenix	Acetone	11185B	1
Hedinger/ Noweda	Acetone	11187B	1
Hedinger	Acetone	11187B	1
Hedinger/Sanacorp	Acetone	11187B	2
Hedinger/Noweda	Acetone	11188B	1
Hedinger	Acetone	11188B	1
Hedinger, Noweda	Acetone	11188B	1
Hedinger	Acetone	11188B	4
Hedinger / Sanacorp	Acetone	11188B	1
Dr.Seeger GmbH.	Acetone	11188B	1
Caelo	Acetone	11193B	1
Hedinger	Acetone	11193B	1
Hedinger/Otto	Acetone	11193B	1
Euro OTC	Acetone	11193B	1
Hedinger	Acetone	11193B	1
Hedinger/Phönix	Acetone	11195B	1
Hedinger	Acetone	11195B	2
PHOENIX	Acetone	11195B	2
Hedinger / Phoenix	Acetone	11195B	1
Seege/ Jenne	Acetone	11199B	2
Hedinger	Acetone	11199B	4
Alliance Health 04.06.2...	Acetone	11199B	1
Hedinger/V.d.L.	Acetone	11199B	1
Hedinger / Noweda	Acetone	11199B	1
Dr.Seeger GmbH	Acetone	11199B	1
Hedinger Phönix	Acetone	11200B	1
Hedinger	Acetone	11200B	1
Hedinger/ Jenne	Acetone	11200B	1
Hedinger/Sanacorp	Acetone	11200B	1
Hedinger / Sanacorp	Acetone	11203B	1
Hed	Acetone	11203B	2
Hedinger / Noweda	Acetone	11203B	1
AHD / Hedinger	Acetone	11203B	1
Hedinger	Acetone	11210B	1
Hedinger	Acetone	11210B	10
Hedinger/Noweda	Acetone	11210B	1
Hedinger / Noweda	Acetone	11210B	2
Sanacorp/Phoenix	Acetone	11210B	1
Hedinger/Gehe	Acetone	11210B	1
Hedinger / Sanacorp	Acetone	11210B	1
Stumpf/Hedinger	Acetone	11210B	1
Apoth.E. Bartel	Acetone	11210B	1
Hedinger	Acetone	11210B	1
Hedinger/Sanacorp	Acetone	11210B	1
AHD 29.10.14 EK 6,70	Acetone	11210B	1
Hedinger / Noweda	Acetone	11215B	1
AHD/18.11.14/Hedinger	Acetone	11215B	1
Hedinger/Phoenix	Acetone	11215B	1
Caelo	Acetone	11215B	1
Hedinger	Acetone	11215B	3
Phönix	Acetone	11215b	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Sanacorp	Acetone	11215B	1
Caelo	Acetone	11218B	1
Alliance health 12.01.15	Acetone	11218B	1
hedinger	Acetone	11218B	2
H.Lamotte/Noweda	Acetone	11210B	2
Caelo	Acetone	11224B	1
Noweda	Acetone	11224B	1
Alliance Healthcare/Hed...	Acetone	11224B	1
Hedinger/Noweda	Acetone	11224B	1
Hedinger/Phoenix	Acetone	11224B	2
Noweda	Acetone	11226B	1
Hedinger/Noweda	Acetone	11226B	2
Hedinger / Noweda	Acetone	11226B	1
Hedinger	Acetone	11226B	1
Hedinger/Kehr	Acetone	11227B	1
Hedinger/Gehe	Acetone	11227B	1
Infectopharm	Acetone	11227B	1
Caelo	Acetone	11227B	1
Sanacorp / Hedinger WE:...	Acetone	11227B	1
Hedinger/Noweda	Acetone	11232B	2
Hedinger	Acetone	11232B	5
Hedinger/Phoenix	Acetone	11232B	1
Hedinger/Sanacorp	Acetone	11234B	1
Phönix	Acetone	11234B	1
Hedinger/Noweda	Acetone	11234B	1
Caelo	Acetone	11234B	2
Hedinger	Acetone	11234B	2
Hedinger/Fiebig	Acetone	11241	1
Hedinger / Noweda	Acetone	11236B	2
Hedinger	Acetone	11241b	1
Phönix	Acetone	11242B	1
Caelo	Acetone	11247B	2
Hedinger	Acetone	11250B	3
Hedinger über Phoenix	Acetone	11247B	1
hedinger	Acetone	11250B	1
Hedinger/Phönix	Acetone	11250B	1
hedinger	Acetone	11252B	2
Caelo	Acetone	11252B	1
Noweda/Hedinger	Acetone	11252B	1
Hedinger	Acetone	11256B	2
Hedinger/Phönix	Acetone	11260B	1
Hedinger	Acetone	11260B	2
Caelo	Acetone	11260B	1
Hedinger/Phoenix	Acetone	11260B	1
HEDINGER/Phönix	Acetone	11260B	1
AHD Hedinger	Acetone	11260B	1
Hedinger	Acetone	112618	1
Hedinger/Noweda	Acetone	11263B	1
Caelo	Acetone	11263B	2
Hedinger	Acetone	112718	1
hedinger	Acetone	11271B	1
Caelo	Acetone	11271B	2
Phoenix	Acetone	11271b	1
Hedinger GmbH	Acetone	11271B	1
Sanacorp	Acetone	11271B	1
Phönix W.E:2.03.18 E...	Acetone	11271B	1
Hedinger	Acetone	11272B	2
Hedinger	Acetone	112788	1
Hedinger/Noweda	Acetone	11272B	1
Hedinger/Sanacorp	Acetone	112788	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Phoenix	Acetone	11272B	1
Hedinger	Acetone	11278B	3
Kehr/A. Hedinger	Acetone	112828	1
Hedinger	Acetone	11282B	2
Caelo	Acetone	11283B	1
hedinger	Acetone	11283B	1
Hedinger/Noweda	Acetone	112868	1
Hedinger	Acetone	11289B	1
Euro OTC	Acetone	11297B	1
Hedinger	Acetone	1261B	1
Hedinger/Alliance	Acetone	112958	1
Fischar	Acetone	2051409	1
Hedinger	Acetone	4061306	1
Fischbar	Acetone	2101205	1
Fischar	Acetone	5041403	1
Fischbar	Acetone	4101302	1
Hedinger	Acetone	8121401	1
Fischar	Acetone	6021303	1
Fischar	Acetone	4091513	1
Hedinger	Acetone	10091304	1
Fischar	Acetone	2051408	1
Fischar/Gehe	Acetone	1016051	1
Otto Fischar GmbH & Co...	Acetone	1015041	1
Fischer/Gehe	Acetone	1015041	1
Otto Fischar	Acetone	1016031	1
Otto Fischar gmbh	Acetone	1017092	1
Fischar / Gehe	Acetone	1015062	1
Fischar	Acetone	1016041	1
Fischbar	Acetone	1016041	1
Fischar	Acetone	1016051	1
Phönix	Acetone	1016111	2
Fischar	Acetone	1017021	1
Fischar	Acetone	1017032	1
fischar / ahd	Acetone	1017092	2
Fischar	Acetone	1017092	1
Fischer/Alliance	Acetone	1018021	1
Hedinger/AHD	Acetone	110958A2010	1
otto fischar	Acetone	12011501	1
Hedinger/ Phoenix	Acetone	1214M1002	1
Hiperscan	Acetone		1
Fischar	Acetone	123456	2
Hedinger	Acetone	16061404	1
Hedinger	Acetone	16081301	1
otto fischar	Acetone	14031502	1
Fischar	Acetone	16101206	1
Hedinger	Acetone	16101402	1
Fischbar	Acetone	17091206	1
Hedinger	Acetone	18011306	1
Hedinger	Acetone	19041301	1
Hedinger	Acetone	220114	1
Fischar	Acetone	22111302	1
Fischbar	Acetone	22091201	1
otto fischar	Acetone	23061503	1
Fischar	Acetone	25101302	1
Hedinger	Acetone	25071411	1
Hedinger	Acetone	2704211002	1
Hedinger	Acetone	2713M1002	1
Fischar	Acetone	28031301	1
	Acetone	281013A	1
Fischar	Acetone	30081301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Gehe	Acetone	311092	1
Fischar	Acetone	4008978161101	1
Fischer/Gehe	Acetone	4201051	1
Kehr	Acetone	4202021	1
Fischar/ Phönix	Acetone	4202022	1
Phönix/Fischer	Acetone	4202052	1
Phönix/Fischar	Acetone	4202052	2
Phönix/Fischar	Acetone	4202061	3
fischer/Anzag	Acetone	4202061	1
Fischar	Acetone	4202061	1
Sanacorp	Acetone	4202062	1
Phönix/Fischar	Acetone	4202062	1
Fischar/Gehe	Acetone	4202101	1
Fischar	Acetone	4202101	1
Fischar/ Phönix	Acetone	4202101	2
Phönix	Acetone	4202112	1
Otto Fischar/Phoenix	Acetone	4202112	1
Fischar	Acetone	4202121	2
Fischer/Sanacorp	Acetone	420301	1
Fischar	Acetone	4203011	4
Fischar/ Phönix	Acetone	4203031	27
Fischar/Gehe	Acetone	4203031	3
fischer/Anzag	Acetone	4203032	2
Fischar	Acetone	4203051	10
Fischar	Acetone	4203032	1
Fischar	Acetone	4203091	1
Otto Fischer/ Gehe	Acetone	4203101	1
Fischar/Noweda	Acetone	4203101	1
Otto Fischar	Acetone	4203101	1
Sanacorp/ Fischar	Acetone	4203111	1
Fischar/Gehe	Acetone	4203101	1
Fischar/ Now	Acetone	4203111	1
Fischar/Phönix	Acetone	4204012	2
Fischar/Gehe	Acetone	4204031	1
Fischar	Acetone	4204031	3
otto fischar/gehe	Acetone	4204032	1
Phönix Gotha	Acetone	4204032	1
Anzag/ Otto Fischar Gmb...	Acetone	4204032	1
Fischar	Acetone	4204032	1
Euro OTC	Acetone	4204032	1
Otto Fischar GmbH&Co.KG...	Acetone	4204032	1
Fischar	Acetone	4204051	1
Otto Fischar GmbH	Acetone	4204032	1
Otto Fischar GmbH & Co...	Acetone	4204051	1
Otto Fischbar	Acetone	4204052	1
Fischar	Acetone	4204052	1
Otto Fischer/ Gehe	Acetone	4204091	1
Fischar	Acetone	4204071	1
Otto Fischar	Acetone	4204091	1
Otto Fischar/Phönix	Acetone	4205012	1
Fischar/Gehe	Acetone	4204071	1
Fischar	Acetone	4204072	1
Fischar	Acetone	4205011	1
fischar	Acetone	4205031	1
Fischar/Phoenix	Acetone	4205031	1
Caelo	Acetone	4205062	1
Fischer	Acetone	4205032	1
Fischar	Acetone	4205061	1
Fischar	Acetone	4205031	1
Fischer/Gehe 3.11.15	Acetone	4205062	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hubertus Apotheke	Acetone	4205062	1
Fischar/Gehe	Acetone	4205071	1
Fischar/Noweda	Acetone	4205062	1
Fischar	Acetone	4205071	1
Fischar/Phönix	Acetone	4205072	1
Otto Fischar GmbH & Co....	Acetone	4205071	1
Otto Fischar GmbH / Gehe	Acetone	4205092	1
Fischar/Gehe	Acetone	4205092	1
Caesar & Loretz GmbH Hi...	Acetone	420511	2
Caesar & Loretz GmbH Hi...	Acetone	4205111	1
Fischar	Acetone	4206011	2
Fischar/Alliance Heathc...	Acetone	4206011	1
Otto Fischar GmbH & Co....	Acetone	4206012	1
Fischar	Acetone	4206012	2
Otto Fischar / Noweda	Acetone	4206012	1
Fischar	Acetone	4206041	1
Fischar/Phönix	Acetone	4206041	1
Ebert KKH 29.09.2016	Acetone	4206042	1
fischer	Acetone	4206061	1
Otto Fischer	Acetone	4206061	1
Sanacorp	Acetone	4206091	1
Fischar GmbH/Krieger	Acetone	4206071	1
Fischar/Noweda	Acetone	4206092	1
Fischar/Noweda	Acetone	4207011	1
Fischar / Gehe	Acetone	4207013	1
Fischar/Phoenix	Acetone	4207013	1
Otto Fischar GmbH	Acetone	4207072	1
Fischar / Sanacorp	Acetone	4207041	1
Fischer/Fiebig	Acetone	4207102	1
fishar/Phoenix	Acetone	4208022	2
Sanacorp	Acetone		2
Fischar	Acetone	45204052	1
Caelo	Acetone	59/AKD1	1
Otto Fischar	Acetone	9011/01	1
Otto Fischar	Acetone	90110/1	1
Fischar	Acetone	90110/1	2
Noweda	Acetone	901101	1
Otto Fischar	Acetone	90205/1	1
Caelo	Acetone		1
Fischar Anzag	Acetone	90205/2	1
Fischer	Acetone	A-040613-3	1
Fischar	Acetone	9C205/1	1
ACM	Acetone	KNC0783J13	2
Fischar / Gehe	Acetone		1
ACM	Acetone	KNC1133J14	1
Hedinger/Phönix	Acetone	11173B	2
Holdermann	Acetone	11158B	1
Hedinger/ Anzag	Acetone	111588	1
Hedinger/Phoenix	Acetone	11187B	1
Dr.Seeger/Jenne	Acetone	11188B	1
Hedinger Jenne	Acetone	11193B	1
Caelo	Acetone	11210B	1
Hedinger/AHD	Acetone	11227B	1
sanacorp, 6,64EUR, 20.04...	Acetone	11256B	1
Hedinger/AHD	Acetone	11282B	1
GEHE	Acetone	4202052	1
Fischar	Acetone	4202062	1
Caelo	Acetone	606001008	1

- 90 881 spectra from 501 *Apo-Ident* customers from a total of 26 830 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Acetone* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Acetone* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	239	0	93 635
Type C	0	437	34	90 881

The substance/substance group *Acetone* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.4895 %)
Type C	100.0000 % (> 99.6586 %)	92.7813 % (> 92.1444 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31575	31575	0.00	145.59
31580	31575	5.56	142.28
31793	31575	3.16	144.50
31796	31575	6.35	139.54

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Acorus oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31159-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Acorus oil; Oleum calami rectific.

Special notes

When selecting the *Acorus oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Acorus oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Acorus oil	12356706	31159	40	1404623
Caelo	Acorus oil	15081901	32180	40	20150320
Caelo	Acorus oil	15081901	32348	40	20150320
Caelo	Acorus oil	15081901	33835	40	20150320

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Acorus oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 384 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 80 spectra from 2 reference samples of the substance/substance group *Acorus oil*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Acorus oil	15081902	32347	40
Caelo	Acorus oil	15081901	33251	40

- 93 794 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 0 spectra from 0 *Apo-Ident* customers from 0 batches from the substance/substance group *Acorus oil*.
- These include spectra of independent samples from 0 batches, of which no spectra have been used for the generation of the database.
- 91 352 spectra from 501 *Apo-Ident* customers from a total of 27 009 batches from a further 334 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Acorus oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Acorus oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	80	0	93 794
Type C	0	0	0	91 352

The substance/substance group *Acorus oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 92.5000 %)
Type C	100.0000 % (> 99.6586 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31159	31159	0.00	34.15
32180	31159	4.86	31.34
32348	31159	5.34	29.82
33835	31159	5.17	33.82

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	All-rac-alpha-tocopheryl acetate
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31020-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

All-rac-alpha-tocopheryl acetate; int-rac-alpha-Tocopherylis acetate; Tocopherylis acetatis

Special notes

When selecting the *All-rac-alpha-tocopheryl acetate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *All-rac-alpha-tocopheryl acetate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	All-rac-alpha-to...	12362307	31020	40	1402688
Caelo	All-rac-alpha-to...	13396304	31371	40	20131211
Euro OTC	All-rac-alpha-to...	1411025-01	31646	40	20141216
Euro OTC	All-rac-alpha-to...	1512007	31817	40	20151216*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *All-rac-alpha-tocopheryl acetate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 468 spectra from 13 reference samples of the substance/substance group *All-rac-alpha-tocopheryl acetate*.
- These include spectra of independent samples from 9 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
zu erfassen	All-rac-alpha-tocopheryl ace...	zu erfassen	30299	30
Caelo	All-rac-alpha-tocopheryl ace...	10159809	30323	40
Euro OTC	All-rac-alpha-tocopheryl ace...	1201009-01	30512	40
Euro OTC	All-rac-alpha-tocopheryl ace...	1411025-01	31646 [†]	18
Euro OTC	All-rac-alpha-tocopheryl ace...	1512007	31817 [†]	20
Euro OTC	All-rac-alpha-tocopheryl ace...	1604060	32286	40
Euro OTC	All-rac-alpha-tocopheryl ace...	1606011	32469	40
Euro OTC	All-rac-alpha-tocopheryl ace...	1704014	33446	40
Caelo	All-rac-alpha-tocopheryl ace...	172397	33688	40
Euro OTC	All-rac-alpha-tocopheryl ace...	1710024	33842	40
Caelo	All-rac-alpha-tocopheryl ace...	172397	33886	40
Euro OTC	All-rac-alpha-tocopheryl ace...	1803012	34096	40
Caelo	All-rac-alpha-tocopheryl ace...	180877	34223	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 93 406 spectra from a total of 2156 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 371 spectra from 154 *Apo-Ident* customers from 120 batches from the substance/substance group *All-rac-alpha-tocopheryl acetate*.
- These include spectra of independent samples from 117 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	All-rac-alpha-tocopheryl ace...	11300005	3
Caelo	All-rac-alpha-tocopheryl ace...	13075202	2
Euro OTC	All-rac-alpha-tocopheryl ace...	1310005-01	2
Caelo	All-rac-alpha-tocopheryl ace...	13201102	1
Caelo	All-rac-alpha-tocopheryl ace...	13217201	14
Caelo	All-rac-alpha-tocopheryl ace...	13217202	4
Caelo	All-rac-alpha-tocopheryl ace...	13217205	12
Caesar & Loretz	All-rac-alpha-tocopheryl ace...	13217205	1
Caelo	All-rac-alpha-tocopheryl ace...	13217207	10
Caelo	All-rac-alpha-tocopheryl ace...	13396304	4
Caelo	All-rac-alpha-tocopheryl ace...	13396305	2
Alliance health 05.01.15	All-rac-alpha-tocopheryl ace...	13396305	1
Caelo	All-rac-alpha-tocopheryl ace...	13396306	3
Caelo	All-rac-alpha-tocopheryl ace...	14011512	1
Euro OTC	All-rac-alpha-tocopheryl ace...	140503201	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1405032-01	3
Caelo	All-rac-alpha-tocopheryl ace...	14076805	8
Caelo	All-rac-alpha-tocopheryl ace...	14076806	6
Caelo	All-rac-alpha-tocopheryl ace...	14076812	8
Caelo	All-rac-alpha-tocopheryl ace...	14076808	2
Caelo	All-rac-alpha-tocopheryl ace...	14076813	3
Euro OTC	All-rac-alpha-tocopheryl ace...	1411025-01	7
Caelo	All-rac-alpha-tocopheryl ace...	14363402	2
Caelo	All-rac-alpha-tocopheryl ace...	14363405	4
Caelo	All-rac-alpha-tocopheryl ace...	14363406	5
Caelo	All-rac-alpha-tocopheryl ace...	14363407	4
Noweda 06.04.2016	All-rac-alpha-tocopheryl ace...	14363406	1
Caelo	All-rac-alpha-tocopheryl ace...	14363408	3
Caelo	All-rac-alpha-tocopheryl ace...	14363410	5
Noweda	All-rac-alpha-tocopheryl ace...	14363410	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1504020	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1504020-01	8
Euro OTC	All-rac-alpha-tocopheryl ace...	1504020-02	1
Caelo	All-rac-alpha-tocopheryl ace...	1504020-02	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1512007-01	5

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	All-rac-alpha-tocopheryl ace...	1512007-01	1
Caelo	All-rac-alpha-tocopheryl ace...	15126801	2
Caelo	All-rac-alpha-tocopheryl ace...	15239002	1
Caelo	All-rac-alpha-tocopheryl ace...	15239001	1
Caelo	All-rac-alpha-tocopheryl ace...	15238808	2
Caelo	All-rac-alpha-tocopheryl ace...	15238806	1
Caelo	All-rac-alpha-tocopheryl ace...	15239003	8
Caelo	All-rac-alpha-tocopheryl ace...	15239009	3
Caelo	All-rac-alpha-tocopheryl ace...	15239010	4
Caelo	All-rac-alpha-tocopheryl ace...	15239011	4
Noweda	All-rac-alpha-tocopheryl ace...	15239011	1
Caelo	All-rac-alpha-tocopheryl ace...	15239012	4
DSM Nutritional Prod./...	All-rac-alpha-tocopheryl ace...	15E18-B07	1
GEHE	All-rac-alpha-tocopheryl ace...	15E18-B07-311818	2
Euro/Fiebig	All-rac-alpha-tocopheryl ace...	1606011-02	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1606011-02	1
Gehe	All-rac-alpha-tocopheryl ace...	16128609	1
Caelo	All-rac-alpha-tocopheryl ace...	16128615	3
phönix	All-rac-alpha-tocopheryl ace...	16128615	1
Caelo	All-rac-alpha-tocopheryl ace...	16128625	2
Euro OTC	All-rac-alpha-tocopheryl ace...	171002401	1
Fagron	All-rac-alpha-tocopheryl ace...	225800-0001	1
Caelo	All-rac-alpha-tocopheryl ace...	26011601	1
Euro OTC	All-rac-alpha-tocopheryl ace...	303462	1
Anzag	All-rac-alpha-tocopheryl ace...	82943389	1
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP13396303	2
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP13396303T	5
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP143634003	1
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP14363403T	7
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP15239006T	2
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP11300004	1
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP12362305T	5
Caelo	All-rac-alpha-tocopheryl ace...	611toc	1
Caelo	All-rac-alpha-tocopheryl ace...	12053802	9
Wepa	All-rac-alpha-tocopheryl ace...	AP15239006T	1
Caelo	All-rac-alpha-tocopheryl ace...	70913Toc	1
Caelo	All-rac-alpha-tocopheryl ace...	12319401	2
Caelo	All-rac-alpha-tocopheryl ace...	12053809	1
Caelo	All-rac-alpha-tocopheryl ace...	12319402	4
Caelo	All-rac-alpha-tocopheryl ace...	12362302	14
Caelo	All-rac-alpha-tocopheryl ace...	12362307	18
Euro OTC	All-rac-alpha-tocopheryl ace...	12362307	2
phönix 30.07.2013	All-rac-alpha-tocopheryl ace...	12362307	1
Sanacorp	All-rac-alpha-tocopheryl ace...	12362307	1
Fagron	All-rac-alpha-tocopheryl ace...	12362307	1
Caelo	All-rac-alpha-tocopheryl ace...	12362404	1
Fagron	All-rac-alpha-tocopheryl ace...	12L14-N06	5
Euro OTC	All-rac-alpha-tocopheryl ace...	1301010-01	3
Sanacorp	All-rac-alpha-tocopheryl ace...	1301010-01	1
Großhandel	All-rac-alpha-tocopheryl ace...	13217201	1
Caelo	All-rac-alpha-tocopheryl ace...	1321702	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1411025	1
Caelo	All-rac-alpha-tocopheryl ace...	16091301	1
Caelo	All-rac-alpha-tocopheryl ace...	19071302	1
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP12362305	1
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP11300004T	5
Caesar & Loretz GmbH Hi...	All-rac-alpha-tocopheryl ace...	11300004	1
Fagron	All-rac-alpha-tocopheryl ace...	14B03-B06-294690	2
Fagron	All-rac-alpha-tocopheryl ace...	14G02-B05-303462	2
Fagron	All-rac-alpha-tocopheryl ace...	14G02-B05-3043462	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz GmbH	All-rac-alpha-tocopheryl ace...	15239011	1
Finteler Apotheke	All-rac-alpha-tocopheryl ace...	15E18-B07-326638	1
Audor Pharma	All-rac-alpha-tocopheryl ace...	AP16128607T	1
Caelo	All-rac-alpha-tocopheryl ace...	12053803	6
Fagron	All-rac-alpha-tocopheryl ace...	12I28-B40	2
Caelo	All-rac-alpha-tocopheryl ace...	3312M-02500	1
Caelo	All-rac-alpha-tocopheryl ace...	1813M-02500	1
Noweda	All-rac-alpha-tocopheryl ace...	912303112	1
Fagron	All-rac-alpha-tocopheryl ace...	UT12020065	1
Sanacorp	All-rac-alpha-tocopheryl ace...	Ap12362305T	1
Caelo	All-rac-alpha-tocopheryl ace...	1411025	1
Fagron	All-rac-alpha-tocopheryl ace...	15e18-b07-311818	1
Fagron	All-rac-alpha-tocopheryl ace...	8120	1
Caelo	All-rac-alpha-tocopheryl ace...	10247405	1
Caelo	All-rac-alpha-tocopheryl ace...	11300001	3
Euro OTC	All-rac-alpha-tocopheryl ace...	11300002	1
Caelo	All-rac-alpha-tocopheryl ace...	11300003	3
Euro OTC	All-rac-alpha-tocopheryl ace...	1208020-01	4
Noweda	All-rac-alpha-tocopheryl ace...	1208020-01	1
ANZAG	All-rac-alpha-tocopheryl ace...	12319401	1
Fagron	All-rac-alpha-tocopheryl ace...	12E16-N06	2
Audor Pharma	All-rac-alpha-tocopheryl ace...	12362305T	2
Fagron	All-rac-alpha-tocopheryl ace...	12G12-No7	1
VDL; 24,15EUR; 26.01.2015	All-rac-alpha-tocopheryl ace...	13217206	1
Gehe	All-rac-alpha-tocopheryl ace...	UT10100297	1
Caesar & Loretz GmbH	All-rac-alpha-tocopheryl ace...	13217202	1
Caelo	All-rac-alpha-tocopheryl ace...	14363401	2
Caelo	All-rac-alpha-tocopheryl ace...	11043007	1
EUro/AHD	All-rac-alpha-tocopheryl ace...	1301010-01	1
Sanacorp	All-rac-alpha-tocopheryl ace...	AP11300004T	1
Caelo	All-rac-alpha-tocopheryl ace...	7123623	1
Fagron	All-rac-alpha-tocopheryl ace...	12/28-B40-288705	1
Spinnrad	All-rac-alpha-tocopheryl ace...	1116303113	1
PHOENIX am 27.03.2014	All-rac-alpha-tocopheryl ace...	1310005-01	1
phönix	All-rac-alpha-tocopheryl ace...	1512007-01	1
Stumpf / Hedinger	All-rac-alpha-tocopheryl ace...	15126801	1
Fagron	All-rac-alpha-tocopheryl ace...	15E18-B07-317323	1
Caelo	All-rac-alpha-tocopheryl ace...	16128610	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1606011-01	1
Spinnrad	All-rac-alpha-tocopheryl ace...	217303111	1
Caelo	All-rac-alpha-tocopheryl ace...	17239704	1
Spinnrad	All-rac-alpha-tocopheryl ace...	915303112	1
Caelo	All-rac-alpha-tocopheryl ace...	11300002	3
Caelo	All-rac-alpha-tocopheryl ace...	11043011	1
Caelo	All-rac-alpha-tocopheryl ace...	113000002	1
Anzag ek: 7.71	All-rac-alpha-tocopheryl ace...	11300002	1
Caesar & Loretz GmbH	All-rac-alpha-tocopheryl ace...	12053802	1
Gehe 16.05.2013	All-rac-alpha-tocopheryl ace...	12053809	1
Caelo	All-rac-alpha-tocopheryl ace...	12319042	1
Erika apotheke	All-rac-alpha-tocopheryl ace...	12319402	1
Caelo	All-rac-alpha-tocopheryl ace...	12319403	1
Gehe	All-rac-alpha-tocopheryl ace...	12L14-N06	1
DSM Nutritional Prod./ ...	All-rac-alpha-tocopheryl ace...	12E16-N06	1
Euro OTC	All-rac-alpha-tocopheryl ace...	1511I-02500	1
Caelo	All-rac-alpha-tocopheryl ace...	Ch12053802	2
Caelo	All-rac-alpha-tocopheryl ace...	Ch11300002	1
Fagron	All-rac-alpha-tocopheryl ace...	UT13110302	1
Caelo	All-rac-alpha-tocopheryl ace...	20231012-4	1

- 90 981 spectra from 501 *Apo-Ident* customers from a total of 26 889 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *All-rac-alpha-tocopheryl acetate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *All-rac-alpha-tocopheryl acetate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	468	0	93 406
Type C	0	366	5	90 981

The substance/substance group *All-rac-alpha-tocopheryl acetate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.7179 %)
Type C	100.0000 % (> 99.6586 %)	98.6523 % (> 97.8437 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31020	31020	0.00	79.11
31371	31817	2.73	77.26
31646	31020	2.91	78.48
31817	31817	0.00	78.36

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Angelica root oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30939-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Angelica root oil; Angelica archangelica oleum; Oleum angelicae (e radice)

Special notes

When selecting the *Angelica root oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Angelica root oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Angelica root oil	11352704	30939	50	1403252
Caelo	Angelica root oil	13157610	31632	40	1511618
Caelo	Angelica root oil	163376	33082	40	20161213*
Caelo	Angelica root oil	172038	33909	40	20170728*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 170 spectra of 4 reference samples from the substance/substance group *Angelica root oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 374 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Angelica root oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Angelica root oil	13157610	31632 [†]	20
Caelo	Angelica root oil	15337207	32929	40
Caelo	Angelica root oil	15337209	33513	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2 spectra from 2 *Apo-Ident* customers from 2 batches from the substance/substance group *Angelica root oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Angelica root oil	15337205	1
Caelo	Angelica root oil	11352706	1

- 91 350 spectra from 501 *Apo-Ident* customers from a total of 27 007 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Angelica root oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Angelica root oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	170	0	35 374
Type B	0	89	11	93 774
Type C	0	2	0	91 350

The substance/substance group *Angelica root oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.4706 %)
Type B	100.0000 % (> 99.9761 %)	89.0000 % (> 86.0000 %)
Type C	100.0000 % (> 99.6630 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30939	30939	0.00	26.98
31632	31632	0.00	20.14
33082	33082	0.00	25.93
33909	33909	0.00	20.52

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Anise oil / Star aniseed oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31158-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Anise oil / Star aniseed oil; Anise oil; Oleum anisi; Oleum anisi stellati; Star anise oil

Special notes

When selecting the *Anise oil / Star aniseed oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Anise oil / Star aniseed oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Star anise oil	13056105	31158	40	1405104
Caelo	Star anise oil	13358202	31404	40	20131105
Euro OTC	Star anise oil	1303017-01	31525	40	20130422

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Anise oil	14230101	31596	40	AR-15-FG-004287-01
Caelo	Anise oil	15148901	31785	40	20150518
Euro OTC	Star anise oil	511026	31810	40	20151214*
Euro OTC	Star anise oil	1609007	32775	40	20160919*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 280 spectra of 7 reference samples from the substance/substance group *Anise oil / Star aniseed oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 7 different batches.
- 35 264 spectra from a total of 833 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 350 spectra from 11 reference samples of the substance/substance group *Anise oil / Star aniseed oil*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Star anise oil	10319003	30205	40
Euro OTC	Star anise oil	1303017-01	31525 [†]	20
Caelo	Anise oil	14230101	31596 [†]	20
Caelo	Anise oil	15148901	31785 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Star anise oil	511026	31810 [†]	20
Caelo	Anise oil	16306802	32930	40
Caelo	Star anise oil	170162	33061	30
Caelo	Star anise oil	15337007	33468	40
Caelo	Anise oil	16306803	33507	40
Caelo	Anise oil	16306803	33509	40
Caelo	Anise oil	171500	33908	40

- 93 524 spectra from a total of 2157 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 144 spectra from 68 *Apo-Ident* customers from 60 batches from the substance/substance group *Anise oil / Star aniseed oil*.
- These include spectra of independent samples from 56 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Anise oil	202170A15428901	1
WE:Sanacorp 14.06.2014. . .	Anise oil	13358202	2
Caelo	Anise oil	14230104	1
Caesar & Loretz	Anise oil	14230104	1
Caelo	Anise oil	15085805	1
Caelo	Anise oil	15119401	4
Caelo	Anise oil	15148904	2
Caelo	Anise oil	15148901	3
Caelo	Anise oil	15428901	2
Caesar & Loretz	Anise oil	15148904	1
Caelo	Anise oil	15428906	3
Caelo	Anise oil	15428903	2
Caesar&Loretz/Noweda	Anise oil	16306802	1
Caelo	Anise oil	1612120A/15428901	1
Caelo	Anise oil	16306802	1
Caelo	Anise oil	16306803	2
Caelo	Anise oil	17150001	5
Caesar & Loretz/ Noweda	Anise oil	17150001	1
Taoasis	Anise oil	42803-127870	1
Apotheke Bauer & Cie. oHG	Anise oil	68952-43-2	1
Caelo	Star anise oil	12322503	3
Caelo	Star anise oil	13358209	3
Caelo	Star anise oil	14195202	6
Caelo	Star anise oil	145230101	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Star anise oil	15337007	1
Caelo	Star anise oil	12322507	5
Euro OTC	Star anise oil	1303017-01	1
EuRho/GEHE	Star anise oil	1303017-01	1
Caelo	Star anise oil	13051404	1
Caelo	Star anise oil	13056105	8
Caelo	Star anise oil	13358208	1
Caelo	Star anise oil	13358202	3
Caelo	Star anise oil	14125501	1
Caelo	Star anise oil	14195208	2
Caelo	Star anise oil	14150301	1
Caelo	Star anise oil	14150303	1
Caelo	Star anise oil	14230102	1
Caelo	Star anise oil	14230104	1
Caelo	Star anise oil	15085805	6
Euro OTC	Star anise oil	1502037-01	2
Caelo	Star anise oil	15085801	1
Caelo	Star anise oil	15337020	4
Caelo	Star anise oil	21102014B	1
Caelo	Star anise oil	22122014C	1
Caelo	Star anise oil	22122014D	1
Caelo	Star anise oil	25082015C	1
Caelo	Star anise oil	27041501	1
Caelo	Star anise oil	29081504	1
Caelo	Star anise oil	31032014C	1
Phönix	Star anise oil	1502037-01	1
Caelo	Star anise oil	15085806	2
Beiersdorf/Noweda	Star anise oil	15337007	2
Caelo	Star anise oil	15337013	5
Caelo	Star anise oil	15337014	1
Caleo/Noweda	Star anise oil	15337020	3
Caesar/Noweda	Star anise oil	411I-01953	1
Caelo	Star anise oil	5021403	1
Caelo	Star anise oil	8012014A	1
Taoasis	Star anise oil	106155	1
Caelo	Star anise oil	17016209	1
Caelo	Star anise oil	11352402	5
Noweda	Star anise oil	11352402	1
Phönix	Star anise oil	11352402	1
Caelo	Star anise oil	11352409	8
Caelo	Star anise oil	1135409	1
Caelo	Star anise oil	1232207	1
Anzag 10.12.12	Star anise oil	1206014-01	1
Caelo	Star anise oil	12322501	2
Klenk	Star anise oil	12322503	1
Gehe 14.06.2013	Star anise oil	12322503	1
Caelo	Star anise oil	12322505	1
Caelo	Star anise oil	13056103	1
Fagron	Star anise oil	12322507	1
Caelo	Star anise oil	13358201	1
Caelo	Star anise oil	21111302	1
Caelo	Star anise oil	28051310	1

- 91 208 spectra from 501 *Apo-Ident* customers from a total of 26 949 batches from a further 332 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Anise oil / Star aniseed oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Anise oil / Star aniseed oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	280	0	35 264
Type B	0	350	0	93 524
Type C	0	143	1	91 208

The substance/substance group *Anise oil / Star aniseed oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.8571 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.2857 %)
Type C	100.0000 % (> 99.6586 %)	99.3056 % (> 97.2222 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31158	31158	0.00	52.42
31404	31810	2.31	55.03
31525	31810	4.03	54.28
31596	31596	0.00	61.46
31785	31596	1.41	59.33
31810	31810	0.00	56.25
32775	32775	0.00	56.76

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Aqueous salicylic acid 50% (non-ionic, hydrophil...
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30520-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)

Special notes

When selecting the *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
apomix	Aqueous salicyli...	21A250913	31119	40	1404661
apomix	Aqueous salicyli...	13B291216	32992	40	20161230
apomix	Aqueous salicyli...	06B090617	33493	40	20170609

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
apomix	Aqueous salicylic acid 50% (...)	09A130911	30362	20
apomix	Aqueous salicylic acid 50% (...)	04A050412	30520	40
apomix	Aqueous salicylic acid 50% (...)	05B160517	33492	40
apomix	Aqueous salicylic acid 50% (...)	11B021117	34035	40

- 93 734 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 168 spectra from 37 *Apo-Ident* customers from 83 batches from the substance/substance group *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)*.
- These include spectra of independent samples from 81 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
apomix	Aqueous salicylic acid 50% (...)	1B20016	1
PKH / Noweda	Aqueous salicylic acid 50% (...)	1B260115	2
Apomix	Aqueous salicylic acid 50% (...)	1B260115	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	1B260115	1
PKH Halle/Sanacorp	Aqueous salicylic acid 50% (...)	1B270114	1
Apomix	Aqueous salicylic acid 50% (...)	1C200116	1
PKH Halle/ Phoenix	Aqueous salicylic acid 50% (...)	1C270114	1
Pharmachem/30	Aqueous salicylic acid 50% (...)	38310314	1
PKH GmbH	Aqueous salicylic acid 50% (...)	2C060314	1
Apomix	Aqueous salicylic acid 50% (...)	2C060314	1
Apomix	Aqueous salicylic acid 50% (...)	3B260315	1
apomix/Phönix	Aqueous salicylic acid 50% (...)	1C200116	1
Apomix/Alliance Healthc...	Aqueous salicylic acid 50% (...)	2C230215	1
apomix/Phönix	Aqueous salicylic acid 50% (...)	2B230215	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	3B260315	1
APOMIX / 30	Aqueous salicylic acid 50% (...)	3B310314	1
apomix	Aqueous salicylic acid 50% (...)	3B260315	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	3C310314	1
Apomix	Aqueous salicylic acid 50% (...)	4B200514	1
Phönix/apomix	Aqueous salicylic acid 50% (...)	4B200514	1
PKH GmbH	Aqueous salicylic acid 50% (...)	4B280415	1
Noweda/apomix	Aqueous salicylic acid 50% (...)	4B280415	1
apomix	Aqueous salicylic acid 50% (...)	4C030516	1
Noweda/PKH GmbH	Aqueous salicylic acid 50% (...)	4C280415	1
apomix	Aqueous salicylic acid 50% (...)	4B280415	1
Henry Lamotte Oils GmbH	Aqueous salicylic acid 50% (...)	5B040615	1
Apomix /Noweda	Aqueous salicylic acid 50% (...)	5B040615	2
apomix	Aqueous salicylic acid 50% (...)	5B040615	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	6B080715	2
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	5B270614	1
Apomix / Noweda	Aqueous salicylic acid 50% (...)	6B080715	2
apomix	Aqueous salicylic acid 50% (...)	6B140616	1
PKH GmbH	Aqueous salicylic acid 50% (...)	6B170714	1
apomix PKH GmbH Halle	Aqueous salicylic acid 50% (...)	6B170714	1
Apomix	Aqueous salicylic acid 50% (...)	6B170714	1
Apomix / GEHE	Aqueous salicylic acid 50% (...)	6B1708714	1
Bombastus	Aqueous salicylic acid 50% (...)	6C170714	1
apomix/Phoenix	Aqueous salicylic acid 50% (...)	6C080715	1
Apomix/Phoenix	Aqueous salicylic acid 50% (...)	6C170714	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	6C170714	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	6C170714	1
Apomix/ Phönix	Aqueous salicylic acid 50% (...)	6C170714	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	7B180814	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	7C280715	1
apomix	Aqueous salicylic acid 50% (...)	8B110816	2
apomix	Aqueous salicylic acid 50% (...)	8B250815	2
Apomix/Alliance Healthc...	Aqueous salicylic acid 50% (...)	8B250815	1
Apomix	Aqueous salicylic acid 50% (...)	8B250815	1
Apomix/Noweda	Aqueous salicylic acid 50% (...)	8B250815	1
Apomix	Aqueous salicylic acid 50% (...)	8C170914	1
Apomix	Aqueous salicylic acid 50% (...)	9B051114	1
Apomix	Aqueous salicylic acid 50% (...)	9B240915	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda/Apomix	Aqueous salicylic acid 50% (...)	9C051114	1
Apomix/Alliance Healthc...	Aqueous salicylic acid 50% (...)	9C240915	1
PKH Halle/ Noweda	Aqueous salicylic acid 50% (...)	10B011214	1
Apomix/Phoenix	Aqueous salicylic acid 50% (...)	10B011214	1
Apomix / 11	Aqueous salicylic acid 50% (...)	10B011214	1
apomix/Anzag	Aqueous salicylic acid 50% (...)	10B011214	1
Apomix/AHD	Aqueous salicylic acid 50% (...)	10B261015	2
apomix	Aqueous salicylic acid 50% (...)	10B261015	1
apomix	Aqueous salicylic acid 50% (...)	10C2601015	1
Apomix/Noweda	Aqueous salicylic acid 50% (...)	12b211215	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	11B021215	1
Caelo	Aqueous salicylic acid 50% (...)	11B021215	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	12B211215	1
apomix/Phoenix	Aqueous salicylic acid 50% (...)	12B211215	1
PKH Halle/ Phoenix	Aqueous salicylic acid 50% (...)	22C161013	1
PKH GmbH	Aqueous salicylic acid 50% (...)	25C041213	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	12B211215	1
Fagron	Aqueous salicylic acid 50% (...)	1B200116	1
Caelo	Aqueous salicylic acid 50% (...)	1B200116	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	1B240117	2
apomix	Aqueous salicylic acid 50% (...)	1B240117	1
Apomix/ Sanacorp	Aqueous salicylic acid 50% (...)	2B230216	1
Caelo	Aqueous salicylic acid 50% (...)	1B240117	1
GEHE	Aqueous salicylic acid 50% (...)	2B050318	1
Apomix	Aqueous salicylic acid 50% (...)	2B270217	2
PKH Halle	Aqueous salicylic acid 50% (...)	2B270217	3
apomix	Aqueous salicylic acid 50% (...)	2B270217	1
Apomix	Aqueous salicylic acid 50% (...)	3B160418	1
Apomix /Noweda	Aqueous salicylic acid 50% (...)	3b230316	1
apomix / Alliance Healt...	Aqueous salicylic acid 50% (...)	3B230316	1
Apomix /Noweda	Aqueous salicylic acid 50% (...)	3b280317	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	4B030516	1
apomix	Aqueous salicylic acid 50% (...)	4B030516	1
PKH Halle/Sanacorp	Aqueous salicylic acid 50% (...)	5B250516	2
Apomix /Noweda	Aqueous salicylic acid 50% (...)	5b250516	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	3C230316	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	4B180417	1
Noweda/PKH GmbH	Aqueous salicylic acid 50% (...)	5B250516	1
GEHE	Aqueous salicylic acid 50% (...)	6B090617	1
Apomix	Aqueous salicylic acid 50% (...)	6B140616	3
apomix/Phönix	Aqueous salicylic acid 50% (...)	6B140616	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	7B060716	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	7B060716	1
Apomix / Noweda	Aqueous salicylic acid 50% (...)	7B300617	1
Apomix	Aqueous salicylic acid 50% (...)	7B300617	1
Beiersdorf/Phönix	Aqueous salicylic acid 50% (...)	7B060716	1
Apomix	Aqueous salicylic acid 50% (...)	8B080817	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	8B080817	1
GEHE	Aqueous salicylic acid 50% (...)	7B300617	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	8B110816	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	8B080817	1
Apomix /Noweda	Aqueous salicylic acid 50% (...)	9b060916	1
Apomix	Aqueous salicylic acid 50% (...)	8B110816	1
Apomix/Gehe	Aqueous salicylic acid 50% (...)	9B060916	1
Apomix	Aqueous salicylic acid 50% (...)	9B070917	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	10B121016	1
Apomix	Aqueous salicylic acid 50% (...)	9B060916	1
Apomix	Aqueous salicylic acid 50% (...)	10B121016	1
GEHE	Aqueous salicylic acid 50% (...)	10B121016	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	10B121016	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Aqueous salicylic acid 50% (...)	10B121017	2
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	10B121017	3
Frago/Noweda	Aqueous salicylic acid 50% (...)	11B021117	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	11B041116	1
PKH Halle/ Phönix	Aqueous salicylic acid 50% (...)	11B041116	1
Apomix/ Sanacorp	Aqueous salicylic acid 50% (...)	11B041116	1
GEHE	Aqueous salicylic acid 50% (...)	11B021117	1
Apomix/Phönix	Aqueous salicylic acid 50% (...)	12B251116	1
GEHE	Aqueous salicylic acid 50% (...)	11B041116	2
JENNE	Aqueous salicylic acid 50% (...)	11B041116	1
Apomix / Noweda	Aqueous salicylic acid 50% (...)	12B251116	1
Apomix/Sanacorp	Aqueous salicylic acid 50% (...)	12C211215	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	13B131217	1
Apomix	Aqueous salicylic acid 50% (...)	13B291216	2
Gehe	Aqueous salicylic acid 50% (...)	13B131217	1
GEHE	Aqueous salicylic acid 50% (...)	13B291216	1
Phoenix (PKH Halle)	Aqueous salicylic acid 50% (...)	3A110213	1
Gehe	Aqueous salicylic acid 50% (...)	4A040313	1
Phoenix/ PKH Halle	Aqueous salicylic acid 50% (...)	5A030413	1
PhR/Sanacorp	Aqueous salicylic acid 50% (...)	4A05/13	1
Gehe	Aqueous salicylic acid 50% (...)	11A300513	1
apomix	Aqueous salicylic acid 50% (...)	6A030413	1
Gehe	Aqueous salicylic acid 50% (...)	9A020812	1
Apomix	Aqueous salicylic acid 50% (...)	9A020812	1
PKH GmbH	Aqueous salicylic acid 50% (...)	12A300513	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	12A300513	1
Apomix	Aqueous salicylic acid 50% (...)	12A300513	1
Apomix	Aqueous salicylic acid 50% (...)	19A160913	1
Apomix / Sanacorp	Aqueous salicylic acid 50% (...)	16A250713	1
PKH Halle / Phoenix	Aqueous salicylic acid 50% (...)	21A250913	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	19A160913	1
Gehe	Aqueous salicylic acid 50% (...)	22B161013	1
Klenk	Aqueous salicylic acid 50% (...)	22B161013	1
Apomix / Phönix	Aqueous salicylic acid 50% (...)	22C161013	1
Apomix / 11	Aqueous salicylic acid 50% (...)	25B041213	1
Apomix	Aqueous salicylic acid 50% (...)	22C161013	1

- 91 184 spectra from 501 *Apo-Ident* customers from a total of 26 926 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	120	20	93 734
Type C	0	164	4	91 184

The substance/substance group *Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	85.7143 % (> 83.5714 %)
Type C	100.0000 % (> 99.6586 %)	97.6190 % (> 95.8333 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31119	31119	0.00	178.40
32992	31119	10.23	182.74
33493	31119	11.55	182.34

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Base cream DAC
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30211-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Base cream DAC; Cremor basalis; Unguentum basale

Special notes

When selecting the *Base cream DAC* substance/substance group, the following information is displayed to the user:

- Advice for the identification of base cream DAC: Improper storage of samples can result in identification problems. Please homogenise the sample and let it adapt to room temperature before the measurement.

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Base cream DAC*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Henry Lamotte	Base cream DAC	5215300	30923	70	20120102
Caelo	Base cream DAC	13139901	31043	40	20130513
Euro OTC	Base cream DAC	1304006-01	31252	40	AR-15-FG-014906-01

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Base cream DAC	14052901	31474	30	AR-15-FG-000163-02
Euro OTC	Base cream DAC	1506025	31709	40	20150618*
Euro OTC	Base cream DAC	1506025	31710	40	20150618*
Euro OTC	Base cream DAC	1506026	31711	40	20150618*
Euro OTC	Base cream DAC	1506026	31712	40	20150618*
Euro OTC	Base cream DAC	1506027	31713	40	20150618*
Euro OTC	Base cream DAC	1506027	31714	40	20150618*
Euro OTC	Base cream DAC	1506028	31715	40	20150618*
Euro OTC	Base cream DAC	1506028	31716	40	20150618*
Caelo	Base cream DAC	154204	31964	40	20151218*
Euro OTC	Base cream DAC	39210516	32345	40	20160623
Caelo	Base cream DAC	17028801	34030	40	AR-18-FG-013988-01

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 620 spectra of 15 reference samples from the substance/substance group *Base cream DAC*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 11 different batches.
- 34 924 spectra from a total of 829 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 5412 spectra from 132 reference samples of the substance/substance group *Base cream DAC*.
- These include spectra of independent samples from 107 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Base cream DAC	10096801	30211	11
Henry Lamotte	Base cream DAC	4758200	30254	41
Henry Lamotte	Base cream DAC	5215300	30744	40
Caelo	Base cream DAC	11297201	30891	60
Henry Lamotte	Base cream DAC	5215300	30923 [†]	10
Euro OTC	Base cream DAC	1301023-01	31005	40
Bombastus	Base cream DAC	276100	31059	40
D'Oroderm	Base cream DAC	1346001	31315	40
Caelo	Base cream DAC	14052901	31474 [†]	10
Euro OTC	Base cream DAC	1506025	31709 [†]	20
Euro OTC	Base cream DAC	1506025	31710 [†]	20
Euro OTC	Base cream DAC	1506026	31711 [†]	20
Euro OTC	Base cream DAC	1506026	31712 [†]	20
Euro OTC	Base cream DAC	1506027	31713 [†]	20
Euro OTC	Base cream DAC	1506027	31714 [†]	20
Euro OTC	Base cream DAC	1506028	31715 [†]	20
Euro OTC	Base cream DAC	1506028	31716 [†]	20
Caelo	Base cream DAC	154198	31963	60
Caelo	Base cream DAC	154204	31964 [†]	20
Caelo	Base cream DAC	154201	31965	60
Caelo	Base cream DAC	153911	31966	60
Caelo	Base cream DAC	153910	31967	60
Caelo	Base cream DAC	153909	31968	60
Caelo	Base cream DAC	153905	31969	60
Caelo	Base cream DAC	153904	31970	60
Caelo	Base cream DAC	153903	31971	60
Caelo	Base cream DAC	153902	31972	60
Caelo	Base cream DAC	153578	31973	60
Caelo	Base cream DAC	153869	31974	60
Caelo	Base cream DAC	153872	31975	60
Caelo	Base cream DAC	153868	31976	60
Caelo	Base cream DAC	153576	31977	60
Caelo	Base cream DAC	160051	31978	60
Caelo	Base cream DAC	160052	31979	60
Caelo	Base cream DAC	160053	31980	60
Caelo	Base cream DAC	160054	31981	60
Caelo	Base cream DAC	160055	31982	60
Caelo	Base cream DAC	161021	32236	40
Caelo	Base cream DAC	161347	32237	40
Caelo	Base cream DAC	161348	32238	40
Caelo	Base cream DAC	161349	32239	40
Caelo	Base cream DAC	161356	32240	40
Caelo	Base cream DAC	161357	32241	40
Caelo	Base cream DAC	161357	32242	40
Caelo	Base cream DAC	161357	32243	40
Caelo	Base cream DAC	161358	32244	40
Caelo	Base cream DAC	161359	32245	40
Caelo	Base cream DAC	161360	32246	40
Caelo	Base cream DAC	161361	32247	40
Caelo	Base cream DAC	161362	32248	40
Caelo	Base cream DAC	161363	32249	40
Caelo	Base cream DAC	161456	32250	40
Caelo	Base cream DAC	161038	32251	40
Caelo	Base cream DAC	160556	32252	40
Caelo	Base cream DAC	160588	32253	40
Caelo	Base cream DAC	160559	32254	40
Caelo	Base cream DAC	161020	32255	40
Caelo	Base cream DAC	161020	32256	40
Caelo	Base cream DAC	161020	32257	40
Caelo	Base cream DAC	161021	32258	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Base cream DAC	161021	32259	40
Caelo	Base cream DAC	161022	32260	40
Caelo	Base cream DAC	161022	32261	40
Caelo	Base cream DAC	161501	32368	40
Caelo	Base cream DAC	161522	32369	40
Caelo	Base cream DAC	161460	32521	40
Caelo	Base cream DAC	162378	32522	40
Caelo	Base cream DAC	161361	32523	40
Caelo	Base cream DAC	161499	32524	40
Caelo	Base cream DAC	161554	32525	40
Caelo	Base cream DAC	162204	32526	40
Caelo	Base cream DAC	162205	32527	40
Caelo	Base cream DAC	162379	32528	40
Caelo	Base cream DAC	162380	32529	40
Caelo	Base cream DAC	162381	32530	40
Caelo	Base cream DAC	162382	32531	40
Caelo	Base cream DAC	162383	32532	40
Caelo	Base cream DAC	162403	32533	40
Caelo	Base cream DAC	162406	32534	40
Caelo	Base cream DAC	162407	32535	40
Caelo	Base cream DAC	162407	32536	40
Caelo	Base cream DAC	162407	32537	40
Caelo	Base cream DAC	162407	32538	40
Caelo	Base cream DAC	162410	32539	40
Caelo	Base cream DAC	162416	32540	40
Caelo	Base cream DAC	162417	32541	40
Caelo	Base cream DAC	162420	32542	40
Caelo	Base cream DAC	162422	32543	40
Caelo	Base cream DAC	162617	32544	40
Caelo	Base cream DAC	162617	32545	40
Caelo	Base cream DAC	162617	32546	40
Caelo	Base cream DAC	161500	32547	40
Caelo	Base cream DAC	161524	32548	40
Caelo	Base cream DAC	161558	32549	40
Caelo	Base cream DAC	171542	33380	40
Caelo	Base cream DAC	171533	33381	40
Caelo	Base cream DAC	171474	33382	40
Caelo	Base cream DAC	171532	33383	40
Caelo	Base cream DAC	171473	33384	40
Caelo	Base cream DAC	171082	33385	40
Caelo	Base cream DAC	171543	33386	40
Caelo	Base cream DAC	171004	33387	40
Caelo	Base cream DAC	171039	33388	40
Caelo	Base cream DAC	171476	33389	40
Caelo	Base cream DAC	171005	33390	40
Caelo	Base cream DAC	171557	33391	40
Caelo	Base cream DAC	171038	33392	40
Caelo	Base cream DAC	171475	33393	40
Caelo	Base cream DAC	171477	33394	40
Caelo	Base cream DAC	171170	33395	40
Caelo	Base cream DAC	171240	33397	40
Caelo	Base cream DAC	171534	33526	40
Caelo	Base cream DAC	171535	33527	40
Caelo	Base cream DAC	171536	33528	40
Caelo	Base cream DAC	171537	33529	40
Caelo	Base cream DAC	171538	33530	40
Caelo	Base cream DAC	171539	33531	40
Caelo	Base cream DAC	171540	33532	40
Caelo	Base cream DAC	171545	33533	40
Caelo	Base cream DAC	171546	33534	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Base cream DAC	171547	33535	40
Caelo	Base cream DAC	171548	33536	40
Caelo	Base cream DAC	171549	33537	40
Caelo	Base cream DAC	171550	33538	40
Caelo	Base cream DAC	171541	33539	40
Caelo	Base cream DAC	171551	33540	40
Caelo	Base cream DAC	171552	33541	40
Caelo	Base cream DAC	171554	33542	40
Caelo	Base cream DAC	171555	33543	40
Caelo	Base cream DAC	171556	33544	40
Caelo	Base cream DAC	171544	33546	40
Caelo	Base cream DAC	172259	33767	40

- 88462 spectra from a total of 2053 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 6280 spectra from 449 *Apo-Ident* customers from 2130 batches from the substance/substance group *Base cream DAC*.
- These include spectra of independent samples from 2108 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Phönix	Base cream DAC	6207	1
Caelo	Base cream DAC	1319901	3
Euro OTC	Base cream DAC	2981016	11
Noweda	Base cream DAC	13347201	1
VDL	Base cream DAC	13380901	1
Caelo	Base cream DAC	13380901	11
Caelo	Base cream DAC	13385301	5
Euro OTC	Base cream DAC	1407018-01	7
Caelo	Base cream DAC	14144101	9
Caelo	Base cream DAC	14090801	10
Caelo	Base cream DAC	14178001	6
Caelo	Base cream DAC	14092001	7
Phoenix Goettin	Base cream DAC	14289601	1
Caelo	Base cream DAC	15167001	13
Caelo	Base cream DAC	15040601	8
Caelo	Base cream DAC	15315401	5
Anzag	Base cream DAC	292258	1
Lamotte/AHD	Base cream DAC	589900	1
Caelo	Base cream DAC	12191901	15

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	12249501	3
Euro OTC	Base cream DAC	1301022-01	6
Noweda	Base cream DAC	131250001	1
Caelo	Base cream DAC	13173501	4
Caelo	Base cream DAC	13139801	5
Caelo	Base cream DAC	13198301	5
Caelo	Base cream DAC	13311401	10
Bombastus	Base cream DAC	281334	2
HENRYLAMOTTE/Gehe	Base cream DAC	5408000	1
Phönix Gotha	Base cream DAC	5489800	1
Lamotte/AHD	Base cream DAC	5489900	5
Caelo	Base cream DAC	72/AKD1	1
Gehe	Base cream DAC	276944	1
Gehe	Base cream DAC	278148	3
Caelo	Base cream DAC	273416352701	1
Caelo	Base cream DAC	273717268701	1
Fagron	Base cream DAC	3987	2
Caelo	Base cream DAC	109647	1
Caelo	Base cream DAC	2020216	2
Ensbona	Base cream DAC	2020216	1
Ensbona/Sanacorp	Base cream DAC	1200216	3
Caelo	Base cream DAC	2121403	1
ENSBONA	Base cream DAC	2020216	1
Noweda	Base cream DAC	2981016	1
Euro OTC	Base cream DAC	298116	1
Euro OTC	Base cream DAC	2991016	7
Euro OTC	Base cream DAC	3001016	4
Noweda	Base cream DAC	3011016	2
Euro OTC	Base cream DAC	300116	1
Caelo	Base cream DAC	3011016	1
Euro OTC	Base cream DAC	3011016	2
Apomix	Base cream DAC	3B250616	1
Apomix/Jenne	Base cream DAC	4D061014	1
Ferd. Eimermacher GmbH ...	Base cream DAC	6671116	3
Caelo	Base cream DAC	6671116	1
gehe	Base cream DAC	5A040414	1
Eimermacher	Base cream DAC	8811114	1
Eimmermacher	Base cream DAC	6681116	1
Fagron	Base cream DAC	6879559	1
Noweda/Ensbona	Base cream DAC	8821114	5
Ensbona	Base cream DAC	8811114	1
Sanacorp/ Ferd. Eimerma...	Base cream DAC	8821114	2
Caelo	Base cream DAC		20
Euro OTC	Base cream DAC	11811116	5
Euro OTC	Base cream DAC	11821116	7
Fagron	Base cream DAC	11831116	1
Euro OTC	Base cream DAC	11831116	1
Euro OTC	Base cream DAC	11841116	1
Caelo	Base cream DAC	1214I01309	1
Eimermacher	Base cream DAC	12170214	2
Beyvers GmbH	Base cream DAC	12170717	1
Caelo	Base cream DAC	126328901	1
Euro OTC	Base cream DAC	1301023-01	7
Gehe	Base cream DAC	1302036-01	1
Caelo	Base cream DAC	13017201	4
Caelo	Base cream DAC	13051503	1
Phoenix	Base cream DAC	1306002-01	1
Euro OTC	Base cream DAC	1306030-01	1
Euro OTC	Base cream DAC	1307003-01	2
Caelo	Base cream DAC	1307003-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Base cream DAC	1307034-01	3
Euro OTC	Base cream DAC	1307044-01	3
Euro OTC	Base cream DAC	1307045-01	1
Euro OTC	Base cream DAC	1307047-01	1
Euro OTC	Base cream DAC	130704701	2
Euro OTC	Base cream DAC	1308010-01	1
Caesar & Loretz GmbH Hi...	Base cream DAC	1308013-01	1
Euro OTC	Base cream DAC	1308013-01	4
Euro OTC	Base cream DAC	1308021-01	3
Euro OTC	Base cream DAC	1308022-01	7
Euro OTC	Base cream DAC	13080224-01	1
Euro OTC	Base cream DAC	1308023-01	3
Euro OTC	Base cream DAC	1308024-01	8
Euro OTC	Base cream DAC	1309006-01	4
Eurho/Noweda	Base cream DAC	1309005-1	2
Euro OTC	Base cream DAC	1309007-01	8
Euro OTC	Base cream DAC	1309008-01	1
Euro OTC	Base cream DAC	1309013-01	7
Euro OTG	Base cream DAC	1309008-01	2
Euro OTC	Base cream DAC	1309015-01	3
Euro OTC	Base cream DAC	131007-01	1
Euro OTC	Base cream DAC	1311004-01	3
Noweda	Base cream DAC	1311004-01	1
Euro OTC	Base cream DAC	1311005-01	5
Euro OTC	Base cream DAC	131100801	1
Caelo	Base cream DAC	131100701	1
Euro OTC	Base cream DAC	1311008-01	2
Caelo	Base cream DAC	1311009-01	1
Euro OTC	Base cream DAC	131101601	1
EuRho/Phoenix	Base cream DAC	1311031-01	1
Euro OTC	Base cream DAC	131103101	1
Euro OTC	Base cream DAC	131200101	1
Euro OTC	Base cream DAC	13111004-01	1
Euro OTC	Base cream DAC	1312004-01	1
Euro/Gehe	Base cream DAC	1312002-01	1
Euro OTC	Base cream DAC	1312003-01	1
EuRho/Phönix	Base cream DAC	1312002-01	1
Euro OTC	Base cream DAC	1312020-01	1
Euro OTC	Base cream DAC	1312021-01	1
Euro OTC	Base cream DAC	1312035-01	6
Euro OTC	Base cream DAC	1312037-01	1
Euro OTC	Base cream DAC	1312036-01	1
EuRho/Phönix	Base cream DAC	1312037-01	1
Noweda	Base cream DAC	1312037-01	1
Euro OTC	Base cream DAC	1312038-01	3
Euro OTC	Base cream DAC	1312039-01	6
Caelo	Base cream DAC	13139901	16
Caelo	Base cream DAC	13158701	1
Beyvers GmbH	Base cream DAC	13170717	1
Caelo	Base cream DAC	13198701	10
Caelo	Base cream DAC	13249701	8
Caelo	Base cream DAC	13311901	5
Caelo	Base cream DAC	13311801	8
Sanacorp	Base cream DAC	13311801	1
Caelo	Base cream DAC	13312001	6
noweda	Base cream DAC	13347101	1
Caelo	Base cream DAC	1334701	1
Caelo	Base cream DAC	13347101	16
Caelo	Base cream DAC	13347201	6
Caelo	Base cream DAC	13347301	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	13377201	11
Caelo	Base cream DAC	13377901	5
Caelo	Base cream DAC	13378001	10
Caelo	Base cream DAC	13378201	7
Caelo	Base cream DAC	13378101	6
Caelo	Base cream DAC	1338091	1
Caelo	Base cream DAC	13381001	14
Caelo	Base cream DAC	13381101	2
Caelo	Base cream DAC	13381201	7
Caelo	Base cream DAC	13384901	10
Caelo	Base cream DAC	13385101	6
Caelo	Base cream DAC	13385201	5
Caelo	Base cream DAC	13386001	7
Caelo	Base cream DAC	13386101	3
Caelo	Base cream DAC	13386201	6
D'Oroderm	Base cream DAC	1346001	1
Caelo	Base cream DAC	1346001	4
Noweda 17.10.2014	Base cream DAC	1346001	2
Fagron	Base cream DAC	13E21-T08-003906	12
Fagron	Base cream DAC	13E21-T09-003910	1
Fagron	Base cream DAC	13E21-T04-003894	3
Caelo	Base cream DAC	13E22-T08-003987	1
Fagron	Base cream DAC	13E22t06003981	1
Fagron	Base cream DAC	13E21-T07-003903	2
Fagron	Base cream DAC	13E23T10004023	1
Fragon	Base cream DAC	13E23-T15-005569	1
Fagron	Base cream DAC	13E23-T16-004035	2
Fargon/Anzag	Base cream DAC	13H14-T01-006203	1
Fagron	Base cream DAC	13H14-T03-006212	1
Fagron	Base cream DAC	13H14-T01-006203	2
Fagron	Base cream DAC	13H14T01006203	1
Fagron	Base cream DAC	13H14-T04-006218	3
Fagron	Base cream DAC	13H14-T06-006224	1
Fagron	Base cream DAC	13H14T080	1
Fagron	Base cream DAC	13H14-T09	1
Fagron	Base cream DAC	13K12-T02-008039	5
Caelo	Base cream DAC	13K12T03008044	1
Fagron	Base cream DAC	13K12-T04-008048	1
Fagron	Base cream DAC	13K12-T04008048	1
Caelo	Base cream DAC	13K12-T06-008054	2
Fagron	Base cream DAC	13K12-T06-0004	1
Caelo	Base cream DAC	13K12-t07-008057	1
Fagron	Base cream DAC	13K12-T07-008057	1
Fagron	Base cream DAC	13K12-T08-008063	1
EuRho	Base cream DAC	13K29-B32-288245	1
Euro OTC	Base cream DAC	140102001	1
Bombastus	Base cream DAC	1401020-01	1
Euro OTC	Base cream DAC	1401020-01	4
Euro OTC	Base cream DAC	1401022-01	5
Euro OTC	Base cream DAC	1401021-01	1
Euro OTC	Base cream DAC	1401023-01	6
Euro OTC	Base cream DAC	1401024-01	8
Caelo	Base cream DAC	1401025-01	3
Euro OTC	Base cream DAC	1401025-01	8
Euro OTC	Base cream DAC	1401026-01	3
EuRho	Base cream DAC	1401026-01	2
EuRho/ Ebert	Base cream DAC	1401026-01	1
Euro OTC	Base cream DAC	140102601	1
Caelo	Base cream DAC	14011515	1
Caelo	Base cream DAC	14035901	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Base cream DAC	14035801	1
Caelo	Base cream DAC	14035801	1
Caelo	Base cream DAC	14036001	2
Caelo	Base cream DAC	14036101	4
Caelo	Base cream DAC	14036201	9
Sanacorp 26.01.15 EK 7,89	Base cream DAC	1403801	1
Caelo	Base cream DAC	140409	1
Caelo	Base cream DAC	14052901	15
Caelo	Base cream DAC	14053001	2
Phoenix	Base cream DAC	14053001	2
Caelo	Base cream DAC	14053101	2
Caelo	Base cream DAC	14053401	4
Caelo	Base cream DAC	14053002	1
Caelo	Base cream DAC	14053301	2
Phönix	Base cream DAC	14053501	1
Caesar & Loretz	Base cream DAC	14053401	1
Caelo	Base cream DAC	14053501	1
Caelo	Base cream DAC	14053601	3
Caelo	Base cream DAC	14053701	5
GEHE 15.05.2014	Base cream DAC	14053701	1
Fagron	Base cream DAC	1405-T07-016800	1
Euro OTC	Base cream DAC	1407019-01	1
Euro OTC	Base cream DAC	1407020-01	1
Euro OTC	Base cream DAC	1407021-01	1
Euro OTC	Base cream DAC	140702201	1
Euro OTC	Base cream DAC	1407022-01	3
Euro OTC	Base cream DAC	1407054-01	3
Euro OTC	Base cream DAC	1407023-01	1
Caelo	Base cream DAC	140705401	1
Euro OTC	Base cream DAC	140705501	1
Euro OTC	Base cream DAC	1407055-01	2
Caelo	Base cream DAC	1408009-01	1
Fagron	Base cream DAC	1408-T08-011710	1
Caelo	Base cream DAC	14090401	8
Sanacorp	Base cream DAC	14090401	1
Caelo	Base cream DAC	1409041	1
Euro OTC	Base cream DAC	1409056-01	1
Caelo	Base cream DAC	14090601	9
Gehe 28.05.2014	Base cream DAC	14090601	1
Phoenix 20.06.2014	Base cream DAC	14090601	1
Raben-Apotheke Mildena	Base cream DAC	14090601	1
Caelo	Base cream DAC	14090701	7
Caelo	Base cream DAC	14091401	3
Caelo	Base cream DAC	14091501	10
Sanacorp	Base cream DAC	14091501	2
Caelo	Base cream DAC	14091601	7
GEHE 24.06.2014	Base cream DAC	14091601	1
Caelo	Base cream DAC	14091701	10
Sanacorp	Base cream DAC	14091701	1
Caelo	Base cream DAC	14091801	9
Caelo	Base cream DAC	14091901	6
Cealo / Noweda	Base cream DAC	14092001	1
Europharm/Otto	Base cream DAC	1411030-01	1
Euro OTC	Base cream DAC	1412002-01	4
Euro OTC	Base cream DAC	1411031-01	1
Ferd. Eimermacher GmbH ...	Base cream DAC	14120715	4
Enzbona/Eimermacher	Base cream DAC	14120715	2
Eimermacher/Sanacorp	Base cream DAC	14120715	4
ensbona/Noweda	Base cream DAC	14120715	1
Enzbona/Eimermacher	Base cream DAC	14120715	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Eimermacher/Noweda	Base cream DAC	14120715	2
Ensbona/Eimermacher	Base cream DAC	14130715	1
Eimermacher/Eimermacher	Base cream DAC	14130715	5
Caelo	Base cream DAC	14130715	1
Eimermacher/Ensbona	Base cream DAC	14130715	1
Ensbona	Base cream DAC	14130715	1
Eimermacher	Base cream DAC	14130817	1
Caelo	Base cream DAC	14143801	6
Caelo	Base cream DAC	14143901	6
Phoenix 13.11.14	Base cream DAC	14143801	1
Raben-Apotheke Mildenau	Base cream DAC	14143901	1
Sanacorp 31.07.2014 EK:...	Base cream DAC	14143901	1
Caelo	Base cream DAC	14144001	4
Sanacorp	Base cream DAC	14144001	1
Noweda	Base cream DAC	14144101	1
Sanacorp	Base cream DAC	14144101	1
Caelo	Base cream DAC	14144201	3
Caelo	Base cream DAC	1415001	1
Caelo	Base cream DAC	14151001	7
Caelo	Base cream DAC	14151101	3
Caelo	Base cream DAC	14151301	2
Caelo	Base cream DAC	14151401	6
Caelo	Base cream DAC	14151501	10
Fagron	Base cream DAC	14151501	1
Caelo	Base cream DAC	14151601	7
Caelo	Base cream DAC	14151801	5
Caelo	Base cream DAC	14151901	3
Schmittmann GmbH	Base cream DAC	14151901	1
Caelo	Base cream DAC	14177501	8
Anzag 19.12.14 3,82	Base cream DAC	14177501	1
Caelo	Base cream DAC	14177701	2
Caelo	Base cream DAC	1417780	1
Caelo	Base cream DAC	14177801	2
Anzag	Base cream DAC	14177801	1
Caelo	Base cream DAC	14177901	5
phoenix 01.10.14	Base cream DAC	14178001	1
Caelo	Base cream DAC	14178101	4
Caelo	Base cream DAC	14178201	13
Caelo	Base cream DAC	14178401	5
Caelo	Base cream DAC	14222801	2
Caelo	Base cream DAC	14223201	11
Euro OTC	Base cream DAC	14223101	1
Caelo	Base cream DAC	14222901	1
Caelo	Base cream DAC	14223001	1
Caelo	Base cream DAC	14227601	3
Caelo	Base cream DAC	14227701	4
Sana 19.09.14 EK 7,89	Base cream DAC	14227801	1
Caelo	Base cream DAC	14227801	8
Sanacorp	Base cream DAC	14227901	1
Caelo	Base cream DAC	14227901	4
Caelo	Base cream DAC	14228001	1
Caelo	Base cream DAC	1423101	1
Euro OTC	Base cream DAC	14232001	1
Caelo	Base cream DAC	14232001	6
Caelo	Base cream DAC	14232101	5
Caelo	Base cream DAC	14232201	10
Allianz Healthcare 20.0...	Base cream DAC	14232201	1
Caelo	Base cream DAC	14232301	4
Sana 17.11.14 EK 7,89	Base cream DAC	14232301	1
Caelo	Base cream DAC	14238101	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	14238201	1
Caelo	Base cream DAC	14238301	5
Caelo	Base cream DAC	14238401	2
Caelo	Base cream DAC	14238501	2
Alliance Healthcare	Base cream DAC	14238501	1
Caelo	Base cream DAC	14244501	4
Caelo	Base cream DAC	14263201	2
Sanacorp	Base cream DAC	14244501	1
Caelo	Base cream DAC	14263301	5
Caelo	Base cream DAC	1426331	1
Caelo	Base cream DAC	14263401	1
Caelo	Base cream DAC	14263601	3
Caelo	Base cream DAC	1428201	1
Caelo	Base cream DAC	14274601	1
Caelo	Base cream DAC	14289601	5
Caelo	Base cream DAC	14289701	6
Caelo	Base cream DAC	14289801	11
Sanacorp	Base cream DAC	14289801	1
Caelo	Base cream DAC	14289901	7
Sanacorp 05.01.15 EK 7,89	Base cream DAC	14290001	1
Caelo	Base cream DAC	14290001	6
Phönix 23.01.15	Base cream DAC	14290001	1
Caelo	Base cream DAC	14290101	2
Caelo	Base cream DAC	142990001	1
Caelo	Base cream DAC	14303101	5
Caelo	Base cream DAC	14303301	3
Caelo	Base cream DAC	14303501	4
Caelo	Base cream DAC	14303401	1
Caelo	Base cream DAC	14303801	8
Sana 11.02.15 EK 8,05	Base cream DAC	14303801	1
Phönix	Base cream DAC	14303901	1
Caelo	Base cream DAC	14303901	13
Caelo	Base cream DAC	14325101	2
Caelo	Base cream DAC	14325201	1
Caelo	Base cream DAC	14325401	1
Gehe 25.02.15	Base cream DAC	14325401	1
Caelo	Base cream DAC	14325501	1
Dr. Wolff	Base cream DAC	14325601	1
Caelo	Base cream DAC	14325801	2
Caelo	Base cream DAC	14325701	1
Caelo	Base cream DAC	14325601	2
Caelo	Base cream DAC	14325901	1
Caelo	Base cream DAC	14326701	12
Caelo	Base cream DAC	14326801	5
Caelo	Base cream DAC	14326901	4
Caelo	Base cream DAC	14327001	6
Caelo	Base cream DAC	14327101	3
Caelo	Base cream DAC	14327301	2
Caelo	Base cream DAC	14327401	3
Caelo	Base cream DAC	14327501	3
Sanacorp	Base cream DAC	14327501	1
Caelo	Base cream DAC	14327701	7
Caelo	Base cream DAC	14327601	3
Caelo	Base cream DAC	14327801	2
Caelo	Base cream DAC	14327901	7
Caelo	Base cream DAC	14328001	2
Caelo	Base cream DAC	14368001	9
Caelo	Base cream DAC	14368101	4
Caelo	Base cream DAC	14368201	4
Caelo	Base cream DAC	14374201	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	14374301	2
Caelo	Base cream DAC	14374401	3
Caelo	Base cream DAC	14374501	3
Fagron	Base cream DAC	14A10-T05-009253	1
Ichthyol / Gehe	Base cream DAC	14A10-T06	1
Fagron	Base cream DAC	14A10-T06	1
Fagron	Base cream DAC	14A10-T06-009259	5
Fagron	Base cream DAC	14A13-T01-009276	1
Fagron	Base cream DAC	14A13-T02-009280	1
Fagron	Base cream DAC	14A13-T04009286	1
Fagron	Base cream DAC	14A13-T04-009286	1
Fagron	Base cream DAC	14A13-T05-009289	1
Fagron	Base cream DAC	14A13-T07-009298	1
Fagron	Base cream DAC	14A13-T06	1
Fagron	Base cream DAC	14B28T04010673	1
Fagron	Base cream DAC	14D085T045011692	1
Fagron	Base cream DAC	14D-08-T01-011683	1
Fagron	Base cream DAC	14D08-T01-011683	3
Phönix	Base cream DAC	14D08-T02-0002	1
Fagron	Base cream DAC	14D08-T02-011686	3
Fagron	Base cream DAC	14d08t02011686	1
Fagron	Base cream DAC	14D08-T03-011689	7
Fagron	Base cream DAC	14D08-T-03-011689	1
Fagron	Base cream DAC	14d08-t04-011632	1
Fagron	Base cream DAC	14D08-T04-011692	8
Fagron	Base cream DAC	14D08-T05-011697	4
Fagron	Base cream DAC	14D08-T05-011698	1
Fagron	Base cream DAC	14D08-T06-011702	2
Fagron	Base cream DAC	14D08-T07-011705	1
Fagron	Base cream DAC	14D08-T08-011710	4
Fagron	Base cream DAC	14D08-T10-011716	2
Noweda	Base cream DAC	14F05-T06-012961	2
Fagron	Base cream DAC	14F05-T06-012961	5
Fagron	Base cream DAC	14F05T06012961	1
Fagron	Base cream DAC	14F05-T07	1
Fagron	Base cream DAC	14F05-T07-012965	4
Caelo	Base cream DAC	14f05-t07-012965	1
Fagron	Base cream DAC	14F05T08	2
Fagron	Base cream DAC	14F05-T22-012975	1
Fagron	Base cream DAC	14F05-T08-012969	2
Fagron	Base cream DAC	14F05-T24	1
Fagron	Base cream DAC	14F05-T24-012983	7
Fagron	Base cream DAC	14F06-T04-013010	1
Fagron	Base cream DAC	14F06-T04	1
Gehe,27.0415	Base cream DAC	14F27-T03-013614	1
Fagron	Base cream DAC	14J06-T08-015886	9
Fagron	Base cream DAC	14I07-T07-015961	1
Fagron	Base cream DAC	14F06-T05-013014	1
Fagron	Base cream DAC	14F06-T04-013010	1
Fagron	Base cream DAC	14I06-T08-015887	1
Fagron	Base cream DAC	14J07-T06-015957	4
Caelo	Base cream DAC	14j07T06015957	1
Fagron	Base cream DAC	14j07-T07-015961	1
Fagron	Base cream DAC	14J07-T07	1
Fagron	Base cream DAC	14J07-T07-015961	2
Fagron	Base cream DAC	14J07-T08-015964	3
Fagron	Base cream DAC	14J07-T09	1
Fagron	Base cream DAC	14J07-T09-015967	3
Fagron	Base cream DAC	14J07-T10-015971	2
Fagron	Base cream DAC	14J07-T10-015972	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	14J07-T11-015975	2
Fagron	Base cream DAC	14J07-T12	1
Fagron	Base cream DAC	14J07-T14-015986	1
Fagron	Base cream DAC	14J07-T13-015981	2
Fagron	Base cream DAC	14J23-T04-016463	2
Fagron	Base cream DAC	14K05T07	1
Fagron	Base cream DAC	14K05-T05-016793	2
Fagron	Base cream DAC	14K05-016793	1
Fagron	Base cream DAC	14K05-T05-061793	1
Fagron	Base cream DAC	14k05-t06-016796	1
Fagron	Base cream DAC	14K05-T06-016796	4
Fagron	Base cream DAC	14K05-T07-016800	7
Fagron	Base cream DAC	14k05t07016800	1
Fagron	Base cream DAC	14k05-t07016800	1
Fagron	Base cream DAC	14K05-T08-016808	1
Sanacorp	Base cream DAC	15000701	1
Caelo	Base cream DAC	15/138/B	1
Caelo	Base cream DAC	15000701	2
Caelo	Base cream DAC	15000801	4
Caelo	Base cream DAC	15000901	6
Caelo	Base cream DAC	15004001	5
Euro OTC	Base cream DAC	1501007-01	1
EuroPharm	Base cream DAC	1501008-01	2
Euro OTC	Base cream DAC	150100801	1
Euro OTC	Base cream DAC	1501008-01	1
Caelo	Base cream DAC	1501019-01	1
Euro OTC	Base cream DAC	1501020-01	1
Euro OTC	Base cream DAC	1501009-01	1
Euro OTC	Base cream DAC	1501019-01	1
Caelo	Base cream DAC	1501039-1	1
Euro OTC	Base cream DAC	1501040-01	1
Europ Ots Pharma	Base cream DAC	1501040-01	1
Euro	Base cream DAC	1501040-01	2
Euro OTC	Base cream DAC	1501050-01	3
Euro OTC	Base cream DAC	1501051-01	4
EuRho	Base cream DAC	1501050-01	1
Euro OTC	Base cream DAC	1501052-01	1
Caelo	Base cream DAC	1501052-01	1
Euro OTC	Base cream DAC	1501055-01	1
Euro OTC	Base cream DAC	1501053-01	2
Euro OTC	Base cream DAC	1501054-01	1
Euro OTC	Base cream DAC	150150-01	1
Fagron	Base cream DAC	15011-T04-020384	1
Euro OTC	Base cream DAC	1502007-01	3
Euro OTC	Base cream DAC	1502008-01	7
EuRho	Base cream DAC	1502008-01	1
Euro OTC	Base cream DAC	1502016	1
Euro OTC	Base cream DAC	1502015-01	4
Euro OTC	Base cream DAC	1502016-01	4
Caelo	Base cream DAC	1502016-01	1
Euro OTC	Base cream DAC	1502017-01	7
EuRho/Phönix	Base cream DAC	1502017-01	1
EuRho/Noweda	Base cream DAC	1502025-01	1
Euro OTC	Base cream DAC	1502025-01	2
Phönix Gotha	Base cream DAC	1502025-01	1
Alliance Healthcare/Hed...	Base cream DAC	1502026-01	2
Euro OTC	Base cream DAC	1502026-01	4
Euro OTC	Base cream DAC	150202601	1
Euro OTC	Base cream DAC	1502027-01	5
Caelo	Base cream DAC	1502027-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Base cream DAC	1502028-01	3
Noweda	Base cream DAC	1502028-01	2
Euro OTC	Base cream DAC	1502029-01	2
EuRho	Base cream DAC	1502029-01	2
Caelo	Base cream DAC	15040701	8
Caelo	Base cream DAC	15040801	3
Ichthyol-Gesellschaft C...	Base cream DAC	15040701	1
Gehe	Base cream DAC	15040801	1
Caelo	Base cream DAC	15040901	4
Caelo	Base cream DAC	15041001	5
Caelo	Base cream DAC	15041101	8
Euro OTC	Base cream DAC	1506026-01	1
Euro OTC	Base cream DAC	1506028-01	4
Euro OTC	Base cream DAC	1506026-011	1
Euro OTC	Base cream DAC	1506045-01	9
Euro OTC	Base cream DAC	150604601	1
Euro OTC	Base cream DAC	1506046-01	2
EuRho / Phoenix	Base cream DAC	1506046-01	1
Euro OTC	Base cream DAC	1506047-01	4
Euro OTC	Base cream DAC	1506053-01	4
Euro OTC	Base cream DAC	1506053	1
Caelo	Base cream DAC	15064301	6
Caelo	Base cream DAC	15064401	6
Caelo	Base cream DAC	15064501	3
Caelo	Base cream DAC	15064601	2
Caelo	Base cream DAC	15064701	11
Caelo	Base cream DAC	150671	1
Sanacorp	Base cream DAC	15070101	1
Caelo	Base cream DAC	15070101	8
Euro OTC	Base cream DAC	1507017-01	1
Euro OTC	Base cream DAC	150701701	1
Euro OTC	Base cream DAC	1507018-01	2
Caelo	Base cream DAC	15070201	8
Euro OTC	Base cream DAC	1508001-01	14
Caelo	Base cream DAC	15071405	1
Caelo	Base cream DAC	15083001	1
Euro OTC	Base cream DAC	1509021-01	14
Caelo	Base cream DAC	15083101	1
Caelo	Base cream DAC	15083201	1
Euro OTC	Base cream DAC	1509021-21	1
Euro OTC	Base cream DAC	1509025-01	1
Euro OTC	Base cream DAC	1509026-01	2
Euro OTC	Base cream DAC	1509028-01	1
Euro OTC	Base cream DAC	1510027-01	1
Basiscreme	Base cream DAC	1510027-01	1
Euro OTC	Base cream DAC	1510030-01	15
Eurho/Noweda	Base cream DAC	1510030-01	2
Euro OTC	Base cream DAC	1511001-01	5
Noweda	Base cream DAC	1511001-01	1
Euro OTC	Base cream DAC	15110011-01	1
Euro OTC	Base cream DAC	1511002-01	11
Fagron	Base cream DAC	1511001-01	1
EuRho	Base cream DAC	1511002-01	2
Euro OTC	Base cream DAC	151100201	1
Bombastus	Base cream DAC	1511002-01	1
Euro OTC	Base cream DAC	1511002-02	2
Euro OTC	Base cream DAC	1512002-01	5
Euro OTC	Base cream DAC	151200201	1
Caelo	Base cream DAC	1513301	1
Caelo	Base cream DAC	15164401	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	151564501	1
Caelo	Base cream DAC	15164301	1
Caelo	Base cream DAC	15164501	2
Caelo	Base cream DAC	15164601	5
Sanacorp	Base cream DAC	15164601	2
Caelo	Base cream DAC	15164701	10
Caelo	Base cream DAC	15166601	4
Caelo	Base cream DAC	15166701	4
Caelo	Base cream DAC	15166801	7
Caelo	Base cream DAC	15166901	3
Caelo	Base cream DAC	15173201	5
Caelo	Base cream DAC	15173301	2
Caelo	Base cream DAC	15173401	2
Caelo	Base cream DAC	15173501	7
Caelo	Base cream DAC	15173601	9
Fagron	Base cream DAC	15173601	1
Caelo	Base cream DAC	15184001	4
Caelo	Base cream DAC	15184101	3
Caelo	Base cream DAC	15184201	1
Caelo	Base cream DAC	15184301	6
Caelo	Base cream DAC	15184301	1
Caelo	Base cream DAC	15184401	14
Caelo	Base cream DAC	15190501	1
Caelo	Base cream DAC	15190601	2
Caelo	Base cream DAC	15190701	2
Caelo	Base cream DAC	15190901	3
Caelo	Base cream DAC	15191001	5
Caelo	Base cream DAC	15191102	5
Caelo	Base cream DAC	15191201	1
Caelo	Base cream DAC	15191202	8
Caelo	Base cream DAC	15191301	2
Caelo	Base cream DAC	15191401	1
Caelo	Base cream DAC	1520601	1
Fagron	Base cream DAC	15241801	1
Caelo	Base cream DAC	15241901	3
Caelo	Base cream DAC	15242001	12
Caelo	Base cream DAC	15250001	3
Caelo	Base cream DAC	15250101	7
Caesar&Loretz,Anzag	Base cream DAC	15250101	1
Sanacorp	Base cream DAC	15250001	1
Caelo	Base cream DAC	15250201	5
Caelo	Base cream DAC	15250301	8
Caelo	Base cream DAC	15250501	13
Gehe	Base cream DAC	15250501	1
Caelo	Base cream DAC	15250601	2
Caelo	Base cream DAC	15250701	2
Caelo	Base cream DAC	15250801	7
Caelo	Base cream DAC	15250902	12
Caelo	Base cream DAC	15272201	15
Caelo	Base cream DAC	15272301	7
Caelo	Base cream DAC	15272402	6
Caesar&Loretz/Caesar&Lo...	Base cream DAC	15272601	1
Caelo	Base cream DAC	15272601	2
Caelo	Base cream DAC	15306101	1
Caelo	Base cream DAC	15306401	9
Caelo	Base cream DAC	15306501	6
Caelo	Base cream DAC	15306701	6
Caelo	Base cream DAC	15306601	2
Caelo	Base cream DAC	15306701	1
Caelo	Base cream DAC	15313501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	15313601	3
Sanacorp	Base cream DAC	15315301	1
Caelo	Base cream DAC	15315301	6
Caelo	Base cream DAC	15318801	5
Caelo	Base cream DAC	15318901	6
Caelo	Base cream DAC	15357601	13
sanacoro	Base cream DAC	15357601	1
Caelo	Base cream DAC	15357701	11
Caelo	Base cream DAC	15357801	6
Caelo	Base cream DAC	15386901	4
Caelo	Base cream DAC	15386801	1
Caelo	Base cream DAC	15387001	3
Bombastus	Base cream DAC	15387001	1
Caelo	Base cream DAC	15387101	6
Sanacorp	Base cream DAC	15387101	1
Caelo	Base cream DAC	15387201	5
Caelo	Base cream DAC	15390401	2
Caelo	Base cream DAC	15390501	3
Caelo	Base cream DAC	15390701	1
Caelo	Base cream DAC	15390601	2
Caelo	Base cream DAC	15390801	6
Sanacorp	Base cream DAC	15390801	2
Caelo	Base cream DAC	15390901	1
Caelo	Base cream DAC	15391001	3
Caelo	Base cream DAC	15391101	2
Caelo	Base cream DAC	1539101	1
Caelo	Base cream DAC	15391201	4
Caelo	Base cream DAC	1541901	1
Caelo	Base cream DAC	15419901	4
Caelo	Base cream DAC	15420001	9
Caelo	Base cream DAC	15420201	10
Caelo	Base cream DAC	15420101	2
WE:16.02.2016	Base cream DAC	15420201	1
Caelo	Base cream DAC	15420401	10
Caelo	Base cream DAC	1542041	1
Euro OTC	Base cream DAC	15502028-01	1
Caelo	Base cream DAC	1587201	1
Fagron	Base cream DAC	1586	1
Fagron	Base cream DAC	15B19-T13-019450	1
Fagron	Base cream DAC	15C11-T04-020383	1
Fagron	Base cream DAC	15C11-T04-020384	3
Fagron	Base cream DAC	15C11T05020388	1
Fagron	Base cream DAC	15C11-T05-020388	7
Bombastus	Base cream DAC	15C-11T05-020388	1
Fagron	Base cream DAC	15C11-T05020388	1
Fagron	Base cream DAC	15C11-T06-020392	2
Fagron	Base cream DAC	15C11-T07	2
Fagron	Base cream DAC	15C11-T07-020396	6
Fagron	Base cream DAC	15C11-T09	1
Fagron	Base cream DAC	15C11-T09-020405	1
Fagron	Base cream DAC	15C11-T11-020411	3
Fagron	Base cream DAC	15c11t07020396	1
Fagron	Base cream DAC	15C11T13020416	1
Fagron	Base cream DAC	15C11-T14-020420	2
Fagron	Base cream DAC	15E27-T02-022881	7
Fagron	Base cream DAC	15C11-To4-020384	1
Fagron	Base cream DAC	15E27-T02-22881	1
Fagron	Base cream DAC	15G09-19-024414	1
Caelo	Base cream DAC	15F30-T14-034792	1
Fagron	Base cream DAC	15G09-T19-024415	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Base cream DAC	15G09-T19-024416	2
Fagron	Base cream DAC	15H04-T05-024975	13
Fagron	Base cream DAC	15H04-T06-024978	6
Fagron	Base cream DAC	15h04-t06-024978	2
Fagron	Base cream DAC	15H04-T07-024983	2
Fagron	Base cream DAC	15H04-T07	1
Fagron	Base cream DAC	15H04T06024978	1
Fagron	Base cream DAC	15H04-T08-024987	2
Fagron	Base cream DAC	15H04-T08	1
Caelo	Base cream DAC	15H04-T08-024987	1
Fagron	Base cream DAC	15H04-T10-024994	2
Fagron	Base cream DAC	15H04-T11-024997	2
Fagron	Base cream DAC	15H04-T12-025000	1
Caelo	Base cream DAC	15H04-T13-025003	2
Fagron	Base cream DAC	15H04-T14-025007	5
Fagron	Base cream DAC	15J13-T03-026895	5
Caelo	Base cream DAC	15J13-T03-026895	1
phönix	Base cream DAC	15j13-t03-026895	1
Fagron	Base cream DAC	15J15-T01-027034	1
Fagron	Base cream DAC	15J13-T04-026901	1
Fagron	Base cream DAC	15117-T01028772	1
Fagron	Base cream DAC	15L17-T01-028772	3
Fagron	Base cream DAC	15117-T03-028775	4
Fagron	Base cream DAC	15L17-T03-028775	7
Fagron	Base cream DAC	15L17-T04	1
Fagron	Base cream DAC	15L17-T04-028778	2
Fagron	Base cream DAC	15L17-T05-028785	2
Caelo	Base cream DAC	16005101	7
Caelo	Base cream DAC	16/078/B	1
Fagron	Base cream DAC	15L17-T06-028791	1
Caelo	Base cream DAC	16005201	3
Bombastus	Base cream DAC	16005201	1
Caelo	Base cream DAC	16005301	6
Caesar&Loretz	Base cream DAC	16005201	1
Frau Brinn-Holtz	Base cream DAC	16005401	1
Caelo	Base cream DAC	16005401	2
Caelo	Base cream DAC	160102801	1
Euro OTC	Base cream DAC	1601027-01	1
Caelo	Base cream DAC	16005501	3
gehe	Base cream DAC	16005501	1
Caelo	Base cream DAC	1601027-01	1
Caelo	Base cream DAC	1601023-01	1
Euro OTC	Base cream DAC	1601055-01	2
Euro OTC	Base cream DAC	1601030-01	1
Euro OTC	Base cream DAC	160123-01	1
Euro OTC	Base cream DAC	1603015-01	2
EuRho/Noweda	Base cream DAC	1603015-01	1
Euro OTC	Base cream DAC	1603016-01	5
Euro OTC	Base cream DAC	160301601	1
Euro OTC	Base cream DAC	1603028-01	9
Euro OTC	Base cream DAC	1603029-01	8
Euro OTC	Base cream DAC	1603030-01	1
Caelo	Base cream DAC	1603037-01	1
Euro OTC	Base cream DAC	1603038-01	1
Caelo	Base cream DAC	1604001-01	1
Euro OTC	Base cream DAC	1603039-01	1
Sanacorp 16.07.2016	Base cream DAC	1604023-01	1
Euro OTC	Base cream DAC	1604040-01	1
Euro OTC	Base cream DAC	1604042-01	2
Euro OTC	Base cream DAC	160404101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	16041001	4
Caelo	Base cream DAC	16041301	3
Caelo	Base cream DAC	16055501	3
Caelo	Base cream DAC	16055601	3
Phoenix	Base cream DAC	16041301	1
Euro OTC	Base cream DAC	1605014-01	1
Caelo	Base cream DAC	16055701	9
Sana 23.03.16 Ek 8,09	Base cream DAC	16055801	1
Caelo	Base cream DAC	16055801	6
Cealo	Base cream DAC	16055801	1
Caelo	Base cream DAC	16055901	3
Caelo	Base cream DAC	1605801	1
Noweda	Base cream DAC	16058801	1
Caelo	Base cream DAC	16058801	12
WE:06.04.2016 Sanacorp	Base cream DAC	16058801	1
Phytocomm	Base cream DAC	16058801	1
Phoenix	Base cream DAC	16059101	1
Caelo	Base cream DAC	16059101	2
Caelo	Base cream DAC	16059201	2
Caelo	Base cream DAC	16061402	1
Caelo	Base cream DAC	160809B	1
Caelo	Base cream DAC	16091901	3
Caelo	Base cream DAC	16101201	7
Caelo	Base cream DAC	16097601	4
Caelo	Base cream DAC	16101301	5
Caelo	Base cream DAC	16101901	7
Caaelo	Base cream DAC	16101901	1
Caelo	Base cream DAC	16102101	1
Caelo	Base cream DAC	16103701	3
Caelo	Base cream DAC	16102201	2
Caelo	Base cream DAC	16103801	6
Caelo	Base cream DAC	16103901	7
Caelo	Base cream DAC	16104001	1
Caelo	Base cream DAC	161111B/16155601	1
Caelo	Base cream DAC	16131501	6
Caelo	Base cream DAC	16128901	3
Caelo	Base cream DAC	16131901	6
Caelo	Base cream DAC	16132001	7
Lamotte	Base cream DAC	16132001	1
Caelo	Base cream DAC	16132101	9
Caelo	Base cream DAC	16132201	2
Caelo	Base cream DAC	16134701	9
Caelo	Base cream DAC	16134901	2
Caelo	Base cream DAC	16135001	10
Caelo	Base cream DAC	16135101	1
Caelo	Base cream DAC	16135701	3
Caelo	Base cream DAC	16135201	2
Caelo	Base cream DAC	16135801	2
Caelo	Base cream DAC	16135901	5
Caelo	Base cream DAC	16136001	3
Caelo	Base cream DAC	16136101	1
Caelo	Base cream DAC	16136201	5
Caelo	Base cream DAC	16136301	6
Caelo	Base cream DAC	16145601	6
Caelo	Base cream DAC	16145701	10
Sanacorp	Base cream DAC	16145701	1
Caelo	Base cream DAC	16145801	7
Caelo	Base cream DAC	16146001	8
Caelo	Base cream DAC	16145901	1
Caelo	Base cream DAC	16149901	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Base cream DAC	1614-T02-031202	1
Caelo	Base cream DAC	16150101	1
Caelo	Base cream DAC	16150201	1
Caeolo/Noweda	Base cream DAC	16152101	1
Caelo	Base cream DAC	16150301	1
Caelo	Base cream DAC	16152101	2
Caelo	Base cream DAC	16152201	9
WE:Sanacorp 10.08.2016...	Base cream DAC	16152201	1
Caelo	Base cream DAC	16152401	3
Caelo	Base cream DAC	161555601	1
WE:Sanacorp26.09.2016...	Base cream DAC	16155401	1
Caelo	Base cream DAC	16155401	3
Caelo	Base cream DAC	16155601	3
Caelo	Base cream DAC	16155801	3
Caelo	Base cream DAC	16155701	1
Caelo	Base cream DAC	16158401	2
Caelo	Base cream DAC	16203301	2
Caelo	Base cream DAC	16205501	3
Gehe	Base cream DAC	16205501	1
Phoenix	Base cream DAC	16205601	1
Caelo	Base cream DAC	16205601	1
Caelo	Base cream DAC	16205801	5
Caelo	Base cream DAC	16206801	4
Caelo	Base cream DAC	16207001	1
Caelo	Base cream DAC	1620901	1
Caelo	Base cream DAC	16214801	7
Caelo	Base cream DAC	16219401	1
Caesar&Loretz	Base cream DAC	16219401	1
Caelo	Base cream DAC	16219501	1
Caelo	Base cream DAC	16220201	3
Caelo	Base cream DAC	16220301	10
Caelo	Base cream DAC	16220401	8
Caelo	Base cream DAC	16220501	2
Caelo	Base cream DAC	16220601	1
Caelo	Base cream DAC	16227001	1
Caelo	Base cream DAC	16227101	6
Caelo	Base cream DAC	16227201	3
Caelo	Base cream DAC	16231901	1
Caelo	Base cream DAC	16232001	5
Ichthyol-Gesellschaft	Base cream DAC	16232001	1
Caeanlo	Base cream DAC	16232001	1
Caelo	Base cream DAC	16232102	3
Bombastus	Base cream DAC	16232201	1
Gehe	Base cream DAC	16232201	1
Caelo	Base cream DAC	16232301	6
WE:Sanacorp 31.10.2016...	Base cream DAC	16232301	1
Caelo	Base cream DAC	16232401	3
Caelo	Base cream DAC	16232601	5
Caelo	Base cream DAC	16232501	1
Caelo	Base cream DAC	16232701	1
Caelo	Base cream DAC	16237901	2
Ichthyol-Gesellschaft	Base cream DAC	16238001	1
Caelo	Base cream DAC	1623801	1
Caelo	Base cream DAC	16238102	2
Caelo	Base cream DAC	16238301	7
Caelo	Base cream DAC	16240301	2
Caelo	Base cream DAC	16240601	3
Caelo	Base cream DAC	16240701	6
Caelo	Base cream DAC	16241001	3
Sanacorp 07.11 Ek.8,04	Base cream DAC	16241001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	16241601	4
Caelo	Base cream DAC	16241701	1
Caelo	Base cream DAC	16241707	1
Caelo	Base cream DAC	16242001	3
Caelo	Base cream DAC	16249801	3
Caelo	Base cream DAC	16242201	2
Caelo	Base cream DAC	16253101	2
Sana 03.12.16 EK 8,04EUR	Base cream DAC	16253101	1
Caelo	Base cream DAC	16253501	2
Caelo	Base cream DAC	16260501	2
Caelo	Base cream DAC	16260701	5
Caelo	Base cream DAC	16260801	5
WESanacorp 17.12.2016 ...	Base cream DAC	16260801	2
Caelo	Base cream DAC	16260901	2
Caelo	Base cream DAC	16261002	3
PHÖNIX	Base cream DAC	16261002	1
Pöckneck Pharma / Gehe	Base cream DAC	16261101	1
Caelo	Base cream DAC	16261101	6
Caelo	Base cream DAC	16261201	2
Caelo	Base cream DAC	16261301	2
Caelo	Base cream DAC	16261401	2
Caelo	Base cream DAC	16261501	3
Caelo/Sanacorp	Base cream DAC	16261501	1
Caelo	Base cream DAC	16261701	2
Caelo	Base cream DAC	16261801	1
Caelo	Base cream DAC	16262001	1
Caelo	Base cream DAC	16261901	1
Caelo	Base cream DAC	16262301	1
Caelo	Base cream DAC	16262101	4
Caelo	Base cream DAC	16318701	5
Caelo	Base cream DAC	16318702	1
Caelo	Base cream DAC	16319301	2
Caelo	Base cream DAC	16319501	2
Caelo	Base cream DAC	16319601	1
Caelo	Base cream DAC	16319401	1
Caelo	Base cream DAC	16325801	6
Caelo	Base cream DAC	16325101	1
Finteler Apotheke	Base cream DAC	16325801	1
Caelo	Base cream DAC	16326001	6
Caelo	Base cream DAC	16325901	1
Caelo	Base cream DAC	16326101	5
Caelo	Base cream DAC	16328901	1
Caelo	Base cream DAC	16329201	3
Caelo	Base cream DAC	16349502	3
Caelo	Base cream DAC	16349501	1
Phoenix	Base cream DAC	16349502	1
Caelo	Base cream DAC	16350101	2
Caelo	Base cream DAC	16350201	4
Caelo	Base cream DAC	16350301	3
Noweda	Base cream DAC	16350201	1
Caesar&Loretz	Base cream DAC	16350501	1
Caelo	Base cream DAC	16352701	6
Sanacorp WE 20.02.17 E...	Base cream DAC	16350301	1
Sanacorp 22.02.17 EK:8,04	Base cream DAC	16352701	1
Caelo	Base cream DAC	16352801	1
Caelo	Base cream DAC	16352901	1
Caelo	Base cream DAC	1650201	1
Fagron	Base cream DAC	16815-T04-030282	1
Caelo	Base cream DAC	1652101	1
Caelo	Base cream DAC	166349502	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Base cream DAC	16A25-T02-029645	1
Fagron	Base cream DAC	16B15-T01-0300274	1
Fagron	Base cream DAC	16b15-f04-030281	1
Fagron	Base cream DAC	16b15-t01-030274	1
Fagron	Base cream DAC	16B15-T01-030274	6
Fagron	Base cream DAC	16b15-T01-030274	1
Fagron	Base cream DAC	16B15-T02-030277	9
Fagron	Base cream DAC	16b15t02-030277	1
Bombastus	Base cream DAC	16B15-T02-030277	1
phönix	Base cream DAC	16b15-t01-030274	1
Caelo	Base cream DAC	16B15-T02-030277	2
Fagron	Base cream DAC	16B15-T04-030281	6
Caelo	Base cream DAC	16B15-T04-030281	4
Fagron	Base cream DAC	16B15-T04-030282	1
phönix	Base cream DAC	16c14-t01-031199	1
Fagron	Base cream DAC	16B15-T05-030285	1
Fagron	Base cream DAC	16C14-T01-031199	1
Fagron	Base cream DAC	16C14-T02-031202	2
Fagron	Base cream DAC	16C14-T02	1
Caelo	Base cream DAC	16C14-T02-031202	1
Fagron	Base cream DAC	16c14t02-031202	1
Fagron	Base cream DAC	16C14-T03	1
Fagron	Base cream DAC	16C14-T03-031205	1
Caelo	Base cream DAC	16C24-T05-031490	2
Fagron	Base cream DAC	16C24-T05-031490	3
Fagron	Base cream DAC	16c24-t05-031490	1
Caelo	Base cream DAC	16c24-T05-031490	1
Fagron GmbH & Co.KG	Base cream DAC	16C24-T05-031490	1
Fagron	Base cream DAC	16c24-t05031490	1
Fagron	Base cream DAC	16C24-T05031490	1
phönix	Base cream DAC	16c24t06031494	1
phönix	Base cream DAC	16C24-T06-031494	1
Fagron	Base cream DAC	16C24-T06-031494	6
Fagron	Base cream DAC	16c24-t06-031494	1
Fagron	Base cream DAC	16D06H07-031930	1
Caelo	Base cream DAC	16C24-To5-031490	1
Fagron	Base cream DAC	16C-T06-031494	1
Fagron	Base cream DAC	16D06-T04-031914	8
Fagron	Base cream DAC	16D06-T04-031914	1
Fagron	Base cream DAC	16D06-T05-031923	5
phönix	Base cream DAC	16D06-T06-031926	1
Fagron	Base cream DAC	16D06-T07-031931	1
Fagron	Base cream DAC	16D06-T06-031926	1
Fagron	Base cream DAC	16D06-T07-0319930	1
Caelo	Base cream DAC	16D06-T10-031951	1
Fagron	Base cream DAC	16D06-T09-031941	1
Fagron	Base cream DAC	16D06-T10-031950	1
Fagron	Base cream DAC	16F30-T13-034789	3
Fagron	Base cream DAC	16F01-T10-033893	1
Fagron	Base cream DAC	16f30-t13-034789	1
Caelo	Base cream DAC	16F30-T13-034789	1
Caelo	Base cream DAC	16F30-T14-034792	2
Fagron	Base cream DAC	16f30-T14-034792	1
Fagron	Base cream DAC	16F30-T14-034792	6
Fagron	Base cream DAC	16F30-T15-034796	1
Fagron	Base cream DAC	16F30-T15-034797	5
phönix	Base cream DAC	16F30-T15-034797	1
Fagron	Base cream DAC	16F30-T16-034800	3
Fagron	Base cream DAC	16G20-T04-035304	1
Fagron	Base cream DAC	16G20-T04-035305	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Base cream DAC	16G20-T05-035308	7
Fagron	Base cream DAC	16G20-T07-035317	1
Fagron	Base cream DAC	16G20-T09-035324	1
Fagron	Base cream DAC	16G20-T08-035321	2
Fagron	Base cream DAC	16g20-t08-035321	1
Fagron	Base cream DAC	16J11-T01-037607	4
Frago/Gehe	Base cream DAC	16J11-T01-037607	1
Fagron	Base cream DAC	16J11-T02-037610	1
Bombastus	Base cream DAC	17000065	2
Fagron	Base cream DAC	16J11-T03-037614	1
phönix	Base cream DAC	16J11-T03-037614	1
Fagron	Base cream DAC	16J11-T04-037619	2
phönix	Base cream DAC	16j11t02037610	1
Bombastus	Base cream DAC	17000069	2
GEHE	Base cream DAC	17000071	1
Bombastus	Base cream DAC	17000072	1
Bombastus	Base cream DAC	17000071	1
Caelo	Base cream DAC	17000102	1
Bombastus	Base cream DAC	17000907	2
Bombastus	Base cream DAC	17000909	1
Bombastus	Base cream DAC	17000914	1
Bombastus	Base cream DAC	17000920	1
Bombastus	Base cream DAC	17000922	1
Caelo	Base cream DAC	17003201	2
Gehe	Base cream DAC	17000922	1
Caelo	Base cream DAC	17003301	1
Caelo	Base cream DAC	17003401	5
Gehe	Base cream DAC	17003501	1
Caelo	Base cream DAC	17003501	3
Caelo	Base cream DAC	17006501	6
Caelo	Base cream DAC	17006601	6
Caelo	Base cream DAC	17006701	2
Bombstus	Base cream DAC	1700907	1
Caelo	Base cream DAC	17009801	1
Caelo	Base cream DAC	17009902	3
Henry Lamotte Oils GMBH	Base cream DAC	170124BC/6003682013	1
Caelo	Base cream DAC	17010201	1
Caelo	Base cream DAC	17014401	3
Caelo	Base cream DAC	17014501	4
Caelo	Base cream DAC	17019501	2
Caelo	Base cream DAC	17014601	1
Caelo	Base cream DAC	170195501	1
Caelo	Base cream DAC	17023301	2
Caelo	Base cream DAC	170216Vi4	1
Caelo	Base cream DAC	17023401	2
Caelo	Base cream DAC	17023501	4
Fagron	Base cream DAC	17023401	1
Caelo	Base cream DAC	17023601	1
Caelo	Base cream DAC	17028801	3
Noweda	Base cream DAC	17028801	1
Caelo	Base cream DAC	17028901	3
Sanacorp 30.03.2017 ek:...	Base cream DAC	17028901	1
Caelo	Base cream DAC	17029001	3
Caeli / Kehr	Base cream DAC	17029001	1
Sanacorp WE 12.04.17 -8...	Base cream DAC	17029001	1
Caelo	Base cream DAC	17029101	2
Caelo	Base cream DAC	17029201	1
Caelo	Base cream DAC	17029301	2
Phoenix	Base cream DAC	17031501	2
Caelo	Base cream DAC	17031601	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Base cream DAC	17031701	2
Caelo	Base cream DAC	17031901	1
Caesar&Loretz	Base cream DAC	17031901	1
Caelo	Base cream DAC	17035701	1
Caelo	Base cream DAC	17035801	4
Henry Lamotte Oils GMBH...	Base cream DAC	170612BC/6004859002	1
Caelo	Base cream DAC	170824CB/17136501	1
Caelo	Base cream DAC	17094401	1
Caelo	Base cream DAC	17090502	1
Caelo	Base cream DAC	17090701	1
Caelo	Base cream DAC	17094501	1
Gehe	Base cream DAC	17090601	1
Caelo	Base cream DAC	170928B/17155101	1
Caelo	Base cream DAC	17097001	5
Caelo	Base cream DAC	17097101	3
Caelo	Base cream DAC	17097202	2
Caelo	Base cream DAC	17097301	3
Gehe	Base cream DAC	17097401	1
Noweda	Base cream DAC	17097401	4
Caelo	Base cream DAC	17097401	4
Caelo	Base cream DAC	17100201	4
Caelo	Base cream DAC	17100301	2
Caelo	Base cream DAC	171003101	1
Caelo	Base cream DAC	171003801	1
Caelo	Base cream DAC	17100401	5
Caelo	Base cream DAC	17100501	3
Fagron	Base cream DAC	17103801	1
Caelo	Base cream DAC	17103701	1
Caelo	Base cream DAC	17103801	1
Caelo	Base cream DAC	17104001	5
Caelo	Base cream DAC	17108101	2
apomix/Noweda	Base cream DAC	17108101	1
Caelo	Base cream DAC	17108201	1
Caelo	Base cream DAC	171114BC/17153802	1
Caelo	Base cream DAC	171128BC/17153802	1
Caelo	Base cream DAC	17116701	12
Caelo	Base cream DAC	17116901	7
Caelo	Base cream DAC	17116801	1
Caelo	Base cream DAC	17121401	3
Caelo	Base cream DAC	17116902	1
Phoenix	Base cream DAC	17121501	1
Caelo	Base cream DAC	17121501	2
Caelo	Base cream DAC	17123901	1
Caelo	Base cream DAC	17124001	4
Caelo	Base cream DAC	17124101	1
Caelo	Base cream DAC	17129901	1
Caelo	Base cream DAC	17133602	1
Phoenix	Base cream DAC	17136301	1
Caelo	Base cream DAC	17136301	3
Caelo	Base cream DAC	17136501	3
Noweda	Base cream DAC	17136601	2
Caelo	Base cream DAC	17147301	3
Caelo	Base cream DAC	17147501	3
Caelo	Base cream DAC	17147701	1
Caelo	Base cream DAC	17153501	1
Caelo	Base cream DAC	17153601	1
Caelo	Base cream DAC	17153701	1
Caesar & Loretz/Noweda	Base cream DAC	17154001	1
Caelo	Base cream DAC	17153901	1
Caelo	Base cream DAC	17154201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	17154401	1
Caelo	Base cream DAC	17154501	1
Caelo	Base cream DAC	17154601	4
Caelo	Base cream DAC	17154801	2
Caelo	Base cream DAC	17154701	1
Caelo	Base cream DAC	17154901	9
Caelo	Base cream DAC	17155001	2
Caesar und Loretz GmbH	Base cream DAC	17155001	1
Caelo	Base cream DAC	17155201	3
Sanacorp WE: 24.08.17 ...	Base cream DAC	17155201	1
Caelo	Base cream DAC	17155701	1
Caelo	Base cream DAC	17155401	1
Caelo	Base cream DAC	17225002	2
Caelo	Base cream DAC	17225201	3
Caelo	Base cream DAC	17225501	1
Phoenix	Base cream DAC	17229201	1
Caelo	Base cream DAC	17225601	2
Sanacop WE:10.11.17 1kg...	Base cream DAC	17225801	1
Caelo	Base cream DAC	17229301	1
Caelo	Base cream DAC	17230101	2
Caelo	Base cream DAC	17229501	1
Caelo	Base cream DAC	17229601	1
Caelo	Base cream DAC	17239901	7
gehe	Base cream DAC	17230001	1
Sanacorp WE:5.01.18	Base cream DAC	17230001	1
Caelo	Base cream DAC	17229801	1
Caelo	Base cream DAC	17240201	1
Caelo	Base cream DAC	17240001	1
Caelo	Base cream DAC	17240301	1
Caelo	Base cream DAC	17240501	1
Caelo	Base cream DAC	17243301	2
Caelo	Base cream DAC	17240801	1
Fagron	Base cream DAC	17243301	1
Phoenix	Base cream DAC	17243501	1
Caelo	Base cream DAC	17243601	5
Caelo	Base cream DAC	17243501	1
Caelo	Base cream DAC	17251501	1
chaelo/gehe	Base cream DAC	17251701	1
Caelo	Base cream DAC	17264201	1
Caelo	Base cream DAC	17264301	2
Caelo	Base cream DAC	17265801	1
Caelo	Base cream DAC	17265901	2
Caelo	Base cream DAC	17264501	1
Caelo	Base cream DAC	17265701	1
Caelo	Base cream DAC	17264601	1
Caelo	Base cream DAC	17266001	3
Sanacorp W.E:17.03.18 ...	Base cream DAC	17266401	1
Caelo	Base cream DAC	17266501	4
Caelo	Base cream DAC	17267001	3
Caelo	Base cream DAC	17266901	1
Sanacoro W.E:01.03.18 ...	Base cream DAC	17266901	1
Caelo	Base cream DAC	17267401	1
Caelo	Base cream DAC	17267501	1
Caelo	Base cream DAC	17267601	2
Caelo	Base cream DAC	17267801	1
Noweda	Base cream DAC	17267801	1
Caelo	Base cream DAC	17268301	2
Caelo	Base cream DAC	17268501	1
Caelo	Base cream DAC	17268601	2
Caelo	Base cream DAC	17268701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	17268801	1
Caelo	Base cream DAC	17268901	1
Caelo	Base cream DAC	17269101	1
Caelo	Base cream DAC	17269201	2
Caelo	Base cream DAC	17349901	1
Caelo	Base cream DAC	17351901	2
Caelo	Base cream DAC	17352101	2
Caelo	Base cream DAC	1768601	1
Caelo	Base cream DAC	17A03-T07-040044	1
Fagron	Base cream DAC	17a03t08040047	1
Caelo	Base cream DAC	17352201	1
Fagron	Base cream DAC	17A03T13040067	1
Fagron	Base cream DAC	17A03-T07-040044	3
Fagron	Base cream DAC	17A03-T08-040047	8
Fagron	Base cream DAC	17A03t10040056	1
Fagron	Base cream DAC	17A03-T10-040056	1
Fagron	Base cream DAC	17A03-T13-040067	3
Frago/Gehe	Base cream DAC	17A03-T11-040060	1
Fagron	Base cream DAC	17A03-T11-040060	1
Fagron	Base cream DAC	17E11-T07-044186	2
Caelo	Base cream DAC	17E11/T08-044210	1
Fagron	Base cream DAC	17F08-T02-045262	4
Fagron	Base cream DAC	17F08T04-045271	1
Phönix	Base cream DAC	17f08-t04-045271	1
Fagron	Base cream DAC		4
Fagron	Base cream DAC	17Fo8-T03-045267	1
Fagron	Base cream DAC	17F08T02045262	1
Fagron	Base cream DAC	17F08-T02-045262	1
Fagron	Base cream DAC	17F08-T03-045267	2
Fagron	Base cream DAC	17H25-T06-047736	5
Caelo	Base cream DAC	17H25-T06-047736	1
Fagron	Base cream DAC	17H25-T07-047741	1
Fagron	Base cream DAC	17K24-T07-050507	2
Fagron	Base cream DAC	17K24-T03-050489	1
Bombastus	Base cream DAC	18000493	1
Bombastus	Base cream DAC	18000957	1
GEHE	Base cream DAC	18000964	1
Bombastus	Base cream DAC	18000972	3
Caelo	Base cream DAC	18001401	1
Caelo	Base cream DAC	18033702	1
Caelo	Base cream DAC	18015101	1
Euro OTC	Base cream DAC	180208BC/86498517	1
Caelo	Base cream DAC	18015301	2
Caesar & Loretz/ Caesar...	Base cream DAC	18037201	1
Caelo	Base cream DAC	18037201	1
Caelo	Base cream DAC	18037402	1
Caelo	Base cream DAC	18055601	2
Caelo	Base cream DAC	18037301	1
Caelo	Base cream DAC	180615010	1
Caelo	Base cream DAC	18091202	1
Caelo	Base cream DAC	18081702	1
Caelo	Base cream DAC	18104701	1
Caelo	Base cream DAC	18105201	2
Caesar & Loretz/ Noweda	Base cream DAC	18105101	1
Caelo	Base cream DAC	18105801	4
Caelo	Base cream DAC	18105101	1
Caelo	Base cream DAC	190814A	1
Caelo	Base cream DAC	18137201	1
Caelo	Base cream DAC	18137801	1
Caelo	Base cream DAC	19081503	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	20011701	1
Caelo	Base cream DAC	19H04T08	1
Caelo	Base cream DAC	20021503	1
Caelo	Base cream DAC	21081407	1
ENSBONA	Base cream DAC	22251214	1
Ensbona/Sanacorp	Base cream DAC	22251214	2
Noweda/Ensbona	Base cream DAC	222512	1
Caelo	Base cream DAC	23061506	1
Caelo	Base cream DAC	23041403	1
Eimermacher	Base cream DAC	23090414	1
Caelo	Base cream DAC	2514U1309	1
Bombastus	Base cream DAC		4
Noweda/Ensbona	Base cream DAC	24180414	4
Caelo	Base cream DAC	26011602	1
Bombastus	Base cream DAC	260416v3	1
Caelo	Base cream DAC	26091405	1
Euro OTC	Base cream DAC	27300917	3
Fagron	Base cream DAC	26F30T14034792	1
Bombastus	Base cream DAC	280078	2
Bombastus	Base cream DAC	278398	5
Bombastus	Base cream DAC	279962	1
Bombastus	Base cream DAC	281243	2
Bombastus	Base cream DAC	282773	2
gehe	Base cream DAC	283915	1
Fagron	Base cream DAC	283915	1
Euro OTC	Base cream DAC	281486	1
Bombastus	Base cream DAC	283918	1
Gehe	Base cream DAC	283974	1
Bombastus	Base cream DAC	283972	2
Bombastus	Base cream DAC	283971	1
Bombastus	Base cream DAC	283973	1
Bombastus	Base cream DAC	283974	2
Bombastus	Base cream DAC	283975	1
Bombastus	Base cream DAC	283976	1
Bombastus	Base cream DAC	285220	1
Bombastus	Base cream DAC	285227	1
Bombastus	Base cream DAC	285228	4
Bombastus	Base cream DAC	286087	5
Bombastus	Base cream DAC	286088	4
Gehe	Base cream DAC	286088	1
Bombastus	Base cream DAC	286089	1
Bombastus	Base cream DAC	286090	2
Bombastus	Base cream DAC	286093	2
Bombastus	Base cream DAC	286094	3
Bombastus	Base cream DAC	286155	2
Fagron	Base cream DAC	286157	1
Bombastus	Base cream DAC	286158	3
Bombastus	Base cream DAC	286159	1
Bombastus	Base cream DAC	286160	3
Bombastus	Base cream DAC	286162	1
Bombastus	Base cream DAC	287110	1
Bombastus	Base cream DAC	287112	2
Bombastus	Base cream DAC	287111	1
Bombastus	Base cream DAC	287114	1
	Base cream DAC	287114	1
Bombastus	Base cream DAC	287115	2
Bombastus	Base cream DAC	287116	2
Bombastus	Base cream DAC	287117	1
Bombastus	Base cream DAC	287196	1
Bombastus	Base cream DAC	287197	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Base cream DAC	287197	1
Bombastus	Base cream DAC	287248	2
Gehe	Base cream DAC	287249	1
Bombastus	Base cream DAC	287251	3
Gehe	Base cream DAC	287868	1
Bombastus	Base cream DAC	287875	2
Bombastus	Base cream DAC	287876	1
Caelo	Base cream DAC	287876	1
Sanacorp	Base cream DAC	287942	1
Bombastus	Base cream DAC	287941	1
Bombastus	Base cream DAC	287944	2
Bombastus	Base cream DAC	288055	4
Bombastus	Base cream DAC	287946	2
Bombastus	Base cream DAC	288054	1
Gehe	Base cream DAC	287947	1
Bombastus	Base cream DAC	288057	1
Bombastus	Base cream DAC	289526	3
	Base cream DAC	289527	1
Bombastus	Base cream DAC	289530	4
Bombastus	Base cream DAC	289531	1
Bombastus	Base cream DAC	289660	1
Bombastus	Base cream DAC	289589	1
Bombastus	Base cream DAC	289663	4
Bombastus	Base cream DAC	289590	1
Sonnenberg-Apotheke	Base cream DAC	289663	1
AHD 11.12.14 EK 8,79	Base cream DAC	289664	1
Bombastus	Base cream DAC	289667	2
Bombastus	Base cream DAC	289665	1
Bombastus	Base cream DAC	289669	1
Fagron	Base cream DAC	289803	1
Bombastus	Base cream DAC	289805	2
Bombastus	Base cream DAC	289802	1
Bombastus	Base cream DAC	289803	1
Bombastus	Base cream DAC	291024	1
Bombastus	Base cream DAC	289804	1
Bombastus	Base cream DAC	290866	1
Gehe	Base cream DAC	291025	1
Bombastus	Base cream DAC	291025	2
Bombastus	Base cream DAC	291026	3
Bombastus	Base cream DAC	291028	1
Bombastus	Base cream DAC	291029	1
Bombastus	Base cream DAC	291032	1
Bombastus	Base cream DAC	291030	1
Phoenix 20.02.15	Base cream DAC	291032	1
Bombastus	Base cream DAC	291034	2
Bombastus	Base cream DAC	291069	1
Bombastus	Base cream DAC	292123	1
Bombastus	Base cream DAC	292122	3
Bombastus	Base cream DAC	292125	4
Bombastus	Base cream DAC	292127	2
Enzbona/Eimermacher	Base cream DAC	29220416	1
Enzbona/Sanacorp	Base cream DAC	29220416	1
Ferd. Eimermacher GmbH ...	Base cream DAC	29220416	3
Caelo	Base cream DAC	29220416	1
Bmbastus/Noweda	Base cream DAC	292256	1
Bombastus	Base cream DAC	292256	3
Bombastus	Base cream DAC	292259	1
Bombastus	Base cream DAC	292260	2
Anzag	Base cream DAC	292260	2
Bombastus	Base cream DAC	292263	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream DAC	292261	1
Bombastus	Base cream DAC	292267	2
Ensbona/Eimermacher	Base cream DAC	29230416	1
Eimermacher/Ensbona	Base cream DAC	29230416	2
Anzag	Base cream DAC	293196	1
Bombastus	Base cream DAC	293200	3
Bombastus	Base cream DAC	293197	1
Bombastus	Base cream DAC	293201	2
Bombastus	Base cream DAC	293196	1
Bombastus	Base cream DAC	293738	1
Bombastus	Base cream DAC	293739	2
Bombastus	Base cream DAC	293795	6
Bombastus	Base cream DAC	293796	1
Bombastus	Base cream DAC	293939	4
Bombastus	Base cream DAC	293800	1
Bombastus	Base cream DAC	293940	1
Bombastus	Base cream DAC	293942	4
Bombastus	Base cream DAC	293943	2
Bombastus	Base cream DAC	293944	3
Bombastus	Base cream DAC	293973	2
Bombastus	Base cream DAC	294829	2
Bombastus	Base cream DAC	294832	2
Bombastus	Base cream DAC	294833	2
Bombastus	Base cream DAC	294831	1
Bombastus	Base cream DAC	294836	1
Bombastus	Base cream DAC	294839	1
Bombastus	Base cream DAC	294838	1
Gehe	Base cream DAC	294834	1
Bombastus	Base cream DAC	294842	3
Bombastus	Base cream DAC	294841	2
Bombastus	Base cream DAC	294843	4
Bombastus	Base cream DAC	295837	1
Bombastus	Base cream DAC	295897	4
Bombastus	Base cream DAC	295899	2
Bombastus	Base cream DAC	295900	1
Bombastus	Base cream DAC	296152	1
Bombastus	Base cream DAC	295901	1
Bombastus	Base cream DAC	296155	5
Bombastus	Base cream DAC	296156	1
AEP	Base cream DAC	296222	1
Bombastus	Base cream DAC	296221	1
Bombastus	Base cream DAC	296224	4
Bombastus	Base cream DAC	296225	6
Bombastus	Base cream DAC	296226	3
Bombastus	Base cream DAC	296227	1
Bombastus	Base cream DAC	296228	1
Bombastus	Base cream DAC	296230	3
Bombastus	Base cream DAC	296229	3
Bombastus	Base cream DAC	296231	3
Bombastus	Base cream DAC	297482	1
Bombastus	Base cream DAC	297484	2
Bombastus	Base cream DAC	297494	4
Bombastus	Base cream DAC	297480	1
Bombastus	Base cream DAC	297496	1
Bombastus	Base cream DAC	297497	1
Bombastus	Base cream DAC	297500	2
Bombastus	Base cream DAC	298170	1
Bombastus	Base cream DAC	298171	1
Bombastus	Base cream DAC	298172	1
Bombastus	Base cream DAC	298173	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream DAC	298256	5
Bombastus	Base cream DAC	298259	3
Caelo	Base cream DAC	298259	2
Bombastus	Base cream DAC	298257	2
Caelo	Base cream DAC	298260	1
Bombastus	Base cream DAC	298260	1
Bombastus	Base cream DAC	298261	2
Bombastus	Base cream DAC	298262	2
Bombastus	Base cream DAC	298263	2
Bombastus	Base cream DAC	298264	2
Bombastus	Base cream DAC	298265	5
Bombastus	Base cream DAC	298266	2
Bombastus	Base cream DAC	298267	1
Euro OTC	Base cream DAC	29870117	2
Euro OTC	Base cream DAC	29860117	2
Gehe	Base cream DAC	29870117	1
Euro OTC	Base cream DAC	29880117	1
Bombastus	Base cream DAC	299319	1
Bombastus	Base cream DAC	299320	3
Bombastus	Base cream DAC	299321	1
Bombastus	Base cream DAC	299322	3
Bombastus	Base cream DAC	299323	4
Bombastus	Base cream DAC	299324	2
Caelo	Base cream DAC	299323	1
Bombastus	Base cream DAC	299945	1
Bombastus	Base cream DAC	299943	1
Bombastus	Base cream DAC	299947	3
Bombastus	Base cream DAC	299950	3
Bombastus	Base cream DAC	299948	3
Caelo	Base cream DAC	2x16261801	1
Bombastus	Base cream DAC	300011	3
Bombastus	Base cream DAC	300013	1
Bombastus	Base cream DAC	300014	5
Bombastus	Base cream DAC	300015	3
Bombastus	Base cream DAC	301148	1
Bombastus	Base cream DAC	301155	3
Caelo	Base cream DAC	301157	1
Bombastus	Base cream DAC	301153	1
Bombastus	Base cream DAC	301157	2
Bombastus	Base cream DAC	301158	3
Bombastus	Base cream DAC	301159	3
Bombastus	Base cream DAC	301160	5
Bombastus	Base cream DAC	301722	1
Bombastus	Base cream DAC	301741	1
Bombastus	Base cream DAC	301743	1
Caelo	Base cream DAC	301743	2
Bombastus	Base cream DAC	301745	6
Bombastus	Base cream DAC	301746	1
Bombastus	Base cream DAC	301754	1
Bombastus	Base cream DAC	301755	1
Bombastus	Base cream DAC	301757	3
Bombastus	Base cream DAC	301756	1
Bombastus	Base cream DAC	301758	3
Bombastus	Base cream DAC	301759	3
Bombastus	Base cream DAC	303055	4
Bombastus	Base cream DAC	303052	2
Bombastus	Base cream DAC	303056	4
Bombastus	Base cream DAC	303058	2
Caelo	Base cream DAC	303058	1
Bombastus	Base cream DAC	303060	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream DAC	303061	2
Caelo	Base cream DAC	303061	1
Bombastus	Base cream DAC	303062	4
Bombastus	Base cream DAC	303063	2
Bombastus	Base cream DAC	303133	1
Bombastus	Base cream DAC	303135	1
Bombastus	Base cream DAC	303414	3
Bombastus	Base cream DAC	303416	1
Bombastus	Base cream DAC	303417	1
Alliance Healthcare	Base cream DAC	303417	1
Bombastus	Base cream DAC	303419	3
Caelo	Base cream DAC	303419	1
Bombastus	Base cream DAC	303422	2
18.10.2016 Phoenix 12,6...	Base cream DAC	303422	1
Bombastus	Base cream DAC	303423	4
Bombastus	Base cream DAC	303420	1
Bombastus	Base cream DAC	303425	1
Bombastus	Base cream DAC	303424	1
Bombastus	Base cream DAC	303428	2
Bombastus	Base cream DAC	303429	5
Bombastus	Base cream DAC	3042362	1
Bombastus	Base cream DAC	304341	1
Bombastus	Base cream DAC	304342	1
Bombastus	Base cream DAC	304343	1
Bombastus	Base cream DAC	304352	3
Bombastus	Base cream DAC	304351	1
Bombastus	Base cream DAC	304353	2
Bombastus	Base cream DAC	304356	1
Bombastus	Base cream DAC	304358	1
Bombastus	Base cream DAC	304354	2
Caelo	Base cream DAC	304352	1
Bombastus	Base cream DAC	304359	2
GEHE	Base cream DAC	304359	1
Bombastus	Base cream DAC	304360	3
Bombastus	Base cream DAC	305710	1
Bombastus	Base cream DAC	305694	1
Bombastus	Base cream DAC	305711	2
Euro OTC	Base cream DAC	305711	1
Bombastus	Base cream DAC	305712	1
GEHE	Base cream DAC	305712	1
Bombastus	Base cream DAC	305714	1
GEHE	Base cream DAC	305713	1
Bombastus	Base cream DAC	305715	3
Bombastus	Base cream DAC	305728	2
Bombastus	Base cream DAC	305730	4
Noweda	Base cream DAC	305715	1
Bombastus	Base cream DAC	305731	2
Bombastus	Base cream DAC	305733	2
Bombastus	Base cream DAC	305734	2
Bombastus	Base cream DAC	305735	1
Bombastus	Base cream DAC	306602	1
Bombastus	Base cream DAC	306611	3
Bombastus	Base cream DAC	306607	1
Bombastus	Base cream DAC	306609	1
Bombastus	Base cream DAC	306613	2
Bombastus	Base cream DAC	306615	4
Bombastus	Base cream DAC	306614	1
Bombastus	Base cream DAC	306617	1
Bombastus	Base cream DAC	306618	2
Caelo	Base cream DAC	306620	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream DAC	306621	3
Bombastus	Base cream DAC	306620	1
Bombastus	Base cream DAC	306623	2
GEHE	Base cream DAC	306621	1
Bombastus	Base cream DAC	306624	2
Bombastus	Base cream DAC	307829	2
Bombastus	Base cream DAC	307830	2
Bombastus	Base cream DAC	307832	1
Bombastus	Base cream DAC	307834	5
Bombastus	Base cream DAC	307963	1
Bombastus	Base cream DAC	307959	1
Bombastus	Base cream DAC	307964	1
GEHE	Base cream DAC	307964	1
Bombastus	Base cream DAC	307965	3
GEHE	Base cream DAC	307965	1
Bombastus	Base cream DAC	307967	2
Bombastus	Base cream DAC	307970	2
Bombastus	Base cream DAC	307968	1
Bombastus	Base cream DAC	308794	1
Bombastus	Base cream DAC	308799	4
Bombastus	Base cream DAC	308795	1
Bombastus	Base cream DAC	308801	2
Bombastus	Base cream DAC	308804	1
Bombastus	Base cream DAC	308805	1
Caelo	Base cream DAC	31031405	1
Caelo	Base cream DAC	31031506	1
GEHE	Base cream DAC	310722	1
GEHE	Base cream DAC	310723	1
Bombastus	Base cream DAC	310723	1
Bombastus	Base cream DAC	310726	3
Bombastus	Base cream DAC	310747	1
Bombastus	Base cream DAC	310746	1
Bombastus	Base cream DAC	310750	1
Bombastus	Base cream DAC	310752	1
Bombastus	Base cream DAC	310751	1
Bombastus	Base cream DAC	310753	5
Bombastus	Base cream DAC	310755	1
Bombastus	Base cream DAC	310756	1
Caelo	Base cream DAC	310756	1
PKH GmbH	Base cream DAC	34567890	1
Euro OTC	Base cream DAC	35980217	1
Euro OTC	Base cream DAC	36920117	1
Euro OTC	Base cream DAC	38540516	12
Caelo	Base cream DAC	39210516	5
Euro OTC	Base cream DAC	3854056	1
Euro OTC	Base cream DAC	39570117	1
Euro OTC	Base cream DAC	395802017	1
Euro OTC	Base cream DAC	39580217	2
Eimermacher/Noweda	Base cream DAC	39590217	1
Euro OTC	Base cream DAC	39600117	1
Euro OTC	Base cream DAC	39620117	3
Ensbona/Sanacorp	Base cream DAC	39590217	1
P.W.Beyvers	Base cream DAC	412/1	2
Beyvers	Base cream DAC	412/1	2
Euro OTC	Base cream DAC	39630217	1
Gehe am 08.04.2014	Base cream DAC	412/1	2
P.W.Beyvers GmbH/Anzag	Base cream DAC	412/2	1
Noweda / Beyvers	Base cream DAC	412/2	1
Fagron	Base cream DAC	412/2	2
Caelo	Base cream DAC	412/2	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Stumpf/Beyvers	Base cream DAC	412/2	1
Beyvers	Base cream DAC	412/3	4
P.W.Beyvers	Base cream DAC	412/3	1
Caelo	Base cream DAC	412/3	1
Noweda/Beyvers	Base cream DAC	412/3	2
Beyvers	Base cream DAC	440/1	3
P.W.Beyvers	Base cream DAC	440/1	1
Bayvers	Base cream DAC	440/1	1
Beyvers	Base cream DAC	440/2	2
Gehe/Beyvers	Base cream DAC	440/2	1
P.W.Beyvers	Base cream DAC	440/2	1
Beyvers	Base cream DAC	440/3	2
Euro OTC	Base cream DAC	44950616	17
Caelo	Base cream DAC	44950616	1
Euro OTC	Base cream DAC	44960616	2
Ensbona/Gehe	Base cream DAC	44960616	1
Euro OTC	Base cream DAC	450380616	1
Euro OTC	Base cream DAC	45330616	16
Ensbona/Eimermacher	Base cream DAC	45150616	1
Euro OTC	Base cream DAC	45340616	8
Caelo	Base cream DAC	45340616	1
Fagron	Base cream DAC	45340616	1
Caelo	Base cream DAC	45348616	1
Caelo	Base cream DAC	45380616	4
Euro OTC	Base cream DAC	45380616	10
Eurho/Noweda	Base cream DAC	45380616	1
Euro OTC	Base cream DAC	46731117	1
Euro OTC	Base cream DAC	49650616	2
Euro OTC	Base cream DAC	49660616	1
Caelo	Base cream DAC	5x17031501	1
Caelo	Base cream DAC	50205/1Gra	1
Noweda	Base cream DAC	50211015	4
Direkt/Phönix	Base cream DAC	50211015	5
Eimermacher	Base cream DAC	50211015	1
Ensbona/ Sanacorp	Base cream DAC	50211015	1
Noweda/ Eimermacher	Base cream DAC	50211015	1
Eimermacher	Base cream DAC	50221015	1
Beyvers / Phönix	Base cream DAC	515/2	2
Phönix Gotha	Base cream DAC	515/1	1
Beyvers	Base cream DAC	515/2	1
P. W. Beyvers GmbH	Base cream DAC	515/2	1
P.W.Beyvers GmbH/ Sanac...	Base cream DAC	515/2	1
AHD/Beyvers	Base cream DAC	515/2	3
Beyvers	Base cream DAC	515/3	1
P.W. Beyvers GmbH	Base cream DAC	515/3	2
Phönix Gotha	Base cream DAC	515/3	1
Bayvers/ AHD	Base cream DAC	515/3	1
Beyvers	Base cream DAC	515/5	1
Euro OTC	Base cream DAC	51560217	2
P.W.Beyvers GmbH / Noweda	Base cream DAC	515/3	1
Noweda	Base cream DAC	51580217	1
Beyvers	Base cream DAC	5153	1
Euro OTC	Base cream DAC	524490616	1
Euro OTC	Base cream DAC	52490616	3
EuRho	Base cream DAC	51580217	1
Euro OTC	Base cream DAC	52500616	1
P.W. Beyvers GmbH / Phö...	Base cream DAC	538/1	2
Phönix Gotha	Base cream DAC	538/1	1
P.W. Beyvers/Noweda	Base cream DAC	538/1	2
Beyvers	Base cream DAC	538/1	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
P.W. Beyvers GmbH	Base cream DAC	538/1	1
AHD/Beyvers	Base cream DAC	538/1	1
Beyvers/AHD	Base cream DAC		1
Beyvers	Base cream DAC	538/2	1
Beyvers/ Phönix	Base cream DAC	538/2	1
Beyvers GmbH	Base cream DAC	538/2	2
Bombastus	Base cream DAC	538/2	2
Beyvers/Noweda	Base cream DAC	538/2	1
Lamotte/ AHDE	Base cream DAC	5407700	1
Henry Lamotte Oils / Gehe	Base cream DAC	5407700	1
Henry Lamotte	Base cream DAC	5489900	2
Beyvers	Base cream DAC	538-1	1
Bayers	Base cream DAC	538/2	1
Lamotte/AHD	Base cream DAC	5490000	2
Ensbona / Ferd. Eimerma...	Base cream DAC	55350716	5
Noweda/ Eimermacher	Base cream DAC	55350716	2
Ensbona	Base cream DAC	55350716	1
Euro OTC	Base cream DAC	55360616	4
Euro OTC	Base cream DAC	553600616	2
EuRho / Phoenix	Base cream DAC	55360616	1
Caelo	Base cream DAC	5537000	8
Phönix	Base cream DAC	55360616	2
Lamotte/Noweda	Base cream DAC	553700	1
ahd/lamotte	Base cream DAC	5537000	2
Henry Lamotte Oils GmbH	Base cream DAC	5537000	12
Fagron	Base cream DAC	5537000	1
Henry Lamotte Oils / Ph...	Base cream DAC	5537000	1
Henry Lamotte/ Noweda	Base cream DAC	5537000	2
Phönix Gotha	Base cream DAC	5537100	1
Henry Lamotte	Base cream DAC	5537100	2
Phoenix 24.03.2014	Base cream DAC	5537200	2
Henry Lamotte/ Noweda	Base cream DAC	5537300	4
Phönix Gotha	Base cream DAC	5537300	1
Lamotte	Base cream DAC	5537300	4
Lamotte/Alliance	Base cream DAC	5537500	2
Henry Lamotte/Anzag	Base cream DAC	5537600	2
Henry Lamotte Oils/Sana...	Base cream DAC	5650600	1
Lamotte/Alliance	Base cream DAC	5537600	1
Euro OTC	Base cream DAC	55850317	1
Phönix Gotha	Base cream DAC	5650600	1
AHD	Base cream DAC	5650700	5
henry lamotte oils gmbh	Base cream DAC	5650900	1
HenryLamotteGmbH/Noweda	Base cream DAC	5650800	1
lamotte/sanacorb	Base cream DAC	5651000	2
Phönix Gotha	Base cream DAC	5650900	1
Henry Lamotte Olis GmbH	Base cream DAC	5651100	1
Lamotte/Anzag	Base cream DAC	5651100	1
Henry Lamotte/Anzag	Base cream DAC	5651200	1
Lamotte/Alliance	Base cream DAC	5651200	1
Phoenix	Base cream DAC	5651300	2
Lamotte/Alliance	Base cream DAC	5651300	1
Noweda/Ensbona	Base cream DAC	56590714	5
Eimermacher	Base cream DAC	57090315	2
Eimermacher	Base cream DAC	57080315	1
Ensbona	Base cream DAC	57090315	4
Ferdinand Eimermacher	Base cream DAC	57090315	3
Ensbona/Sanacorp	Base cream DAC	57090315	3
Noweda/Ensbona	Base cream DAC	57090315	1
Phönix Gotha	Base cream DAC	5787900	1
lamotte/sanacorb	Base cream DAC	5787800	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	5788000	1
Phönix Gotha	Base cream DAC	5788000	1
Lamotte//AHD	Base cream DAC	5787900	1
Caelo	Base cream DAC	5788100	1
Henry Lamotte Oils / No...	Base cream DAC	5788200	1
Henry Lamotte/Anzag	Base cream DAC	5788600	1
Henry Lamotte Oils/Phönix	Base cream DAC	5788700	2
Henry Lamotte/Anzag	Base cream DAC	5788700	3
Lamotte/Sanacorp	Base cream DAC	5788700	2
Caelo	Base cream DAC	5909200	2
Henry Lamotte Oils	Base cream DAC	5909300	2
Sanacorp/Henry Lamotte	Base cream DAC	5909300	1
Caelo	Base cream DAC	5909300	1
Caelo	Base cream DAC	5909500	1
Henry Lamotte	Base cream DAC	5909500	1
Henry Lamotte/Anzag	Base cream DAC	5909500	1
Henry Lamotte/ Noweda	Base cream DAC	5909300	1
Henry Lamotte Oils GmbH	Base cream DAC	5909400	1
Henry Lamotte /Noveda	Base cream DAC	5909600	1
Henry Lamotte/AHD	Base cream DAC	5909700	2
Henry Lamotte Oils GmbH	Base cream DAC	5909600	2
Henry Lamotte Oils / No...	Base cream DAC	5909500	1
Henry Lamotte/Anzag	Base cream DAC	5909800	1
Henry Lamotte	Base cream DAC	5909600	1
Lamotte/Phönix	Base cream DAC	5909600	1
Anzag/Lamotte	Base cream DAC	591000	1
Henry Lamotte Oils/ San...	Base cream DAC	5909700	1
Lamotte/AHD	Base cream DAC	5909600	1
Lamotte	Base cream DAC	590990	1
Henry Lamotte	Base cream DAC	6002237001	2
AHD Lamotte GmbH	Base cream DAC	5910000	1
Henry Lamotte/AHD	Base cream DAC	5909900	1
Henry Lamotte	Base cream DAC	6002237002	2
Lamotte / AHD	Base cream DAC	6002237001	1
Caelo	Base cream DAC	6002237002	1
Lamotte/Gehe	Base cream DAC	6002237002	1
Henry Lamotte Oils GmbH	Base cream DAC	6002237002	1
Lamotte /Noweda	Base cream DAC	6002237002	1
Sana/Henry Lamotte	Base cream DAC	6002237003	1
Lamotte//AHD	Base cream DAC	6002237004	1
Sanacorp/Henry Lamotte	Base cream DAC	6002237005	2
Henry Lamotte Oils / No...	Base cream DAC	6002237005	1
Lamotte/ANZAG	Base cream DAC	6002237008	1
Henry Lamotte Oils	Base cream DAC	6002237005	2
Lamotte/AHD	Base cream DAC	6002237006	1
Lamotte/Sanacorp	Base cream DAC	6002237009	1
Henry Lamotte/Anzag	Base cream DAC	6002237009	2
Lamotte/Sanacorp	Base cream DAC	6002237010	1
Henry Lamotte Oils	Base cream DAC	6002237006	1
Lamotte/AHD	Base cream DAC	600237004	1
Henry/Sanacorp	Base cream DAC	600223705	1
Anzag/Lamotte	Base cream DAC	6003682004	1
Lamotte/Sanacorp	Base cream DAC	6003682005	1
Caelo	Base cream DAC	6003682004	1
Lamotte/ Anzag	Base cream DAC	6003682005	1
Caelo	Base cream DAC	6003682012	1
Lamotte/AHD	Base cream DAC	6003682013	1
Lamotte/AHD	Base cream DAC	6003682012	1
henry lamotte	Base cream DAC	6003682013	2
Caelo	Base cream DAC	6003682015	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte/OG	Base cream DAC	6003682013	1
Caelo	Base cream DAC	6003682014	1
Lamotte/AHD	Base cream DAC	6003682016	1
AH Henry Lamotte GmbH	Base cream DAC	6003682017	2
Alliance/Henry	Base cream DAC	6004859002	2
henry lamotte	Base cream DAC	6004859001	2
Lamotte / Noweda	Base cream DAC	6004859003	1
Henry Lamotte/Noewda	Base cream DAC	6004859001	1
Lamotte/AHD	Base cream DAC	6004859003	1
Lamotte / Noweda	Base cream DAC	6004859002	1
Henry Lamotte/AHD	Base cream DAC	6004859003	2
Lamotte / Noweda	Base cream DAC	6004859004	1
Lamotte/AHD	Base cream DAC	6004859004	1
Gehe	Base cream DAC	6004859005	2
Caelo	Base cream DAC	6004860001	1
henry lamotte	Base cream DAC	6004859005	1
Lamotte / Noweda	Base cream DAC	6004859006	1
AH Henry Lamotte GmbH	Base cream DAC	6004860001	1
Bombastus	Base cream DAC	60090502	2
AHD Lamotte GmbH	Base cream DAC	6004860002	1
Henry Lamotte Oils	Base cream DAC	60260	1
Bombastus	Base cream DAC	60090504	2
Ferdinand Eimermacher	Base cream DAC	61361217	4
Eimermacher/Ensbona	Base cream DAC	63340317	1
Ensbona/Nordwalde	Base cream DAC	63350317	1
Ensbona	Base cream DAC	63350317	2
Ensbona / Sanacorp	Base cream DAC	63340317	2
Sanacorp/Beyvers	Base cream DAC	640/1	3
Ensbona/Sanacorp	Base cream DAC	63350317	1
Beyvers GmbH	Base cream DAC	640/1	1
P.W.Beyvers	Base cream DAC	640/1	2
Beyvers/Noweda	Base cream DAC	640/2	3
Beyvers	Base cream DAC	640/2	1
pharmachem	Base cream DAC	640/2	2
AHD/ Beyvers	Base cream DAC	640/2	2
P.W. Beyvers GmbH/Gehe	Base cream DAC	640/2	1
AHD/ Beyvers	Base cream DAC	640/3	3
Beyvers/Noweda	Base cream DAC	640/3	2
p w Beyvers	Base cream DAC	640/3	1
medphano	Base cream DAC	6500115	3
Euro OTC	Base cream DAC	69500118	1
Euro OTC	Base cream DAC	69520118	3
Euro OTC	Base cream DAC	701228	1
Fagron	Base cream DAC	701228-0002	1
Sanacorp	Base cream DAC	72/AKD2	1
Ensbona/Sanacorp	Base cream DAC	72240417	2
Bayvers	Base cream DAC	731/1	3
Ensbona/phönix	Base cream DAC	72240417	1
AHD/ Beyvers	Base cream DAC	731/1	1
Sanacorp/Beyvers	Base cream DAC	731/2	3
Bayvers	Base cream DAC	731/3	4
Sanacorp/Beyvers	Base cream DAC	731/3	2
P. W. Beyvers GmbH	Base cream DAC	731/4	1
beyvers	Base cream DAC	731/4	1
Caelo	Base cream DAC	74/AKD2	1
Sanacorp	Base cream DAC	75/AKD4	1
Sanacorp	Base cream DAC	75/AKD6	1
Caelo	Base cream DAC	75/AKD7	1
Ensbona	Base cream DAC	79000916	3
Sanacorp	Base cream DAC	78/AKD12	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	79/AKD16	1
AHD Lamotte GmbH	Base cream DAC	800287008	1
henry lamotte	Base cream DAC	8003325003	2
Lamotte / Noweda	Base cream DAC	8003325003	1
Lamotte / Noweda	Base cream DAC	8003325004	2
Euro OTC	Base cream DAC	8003325004	1
henry lamotte	Base cream DAC	8003325006	1
Lamotte / Noweda	Base cream DAC	8003325007	1
AHD Lamotte GmbH	Base cream DAC	8003325009	1
Ferdinand Eimermacher	Base cream DAC	80340218	2
Henry/Sanacorp	Base cream DAC	81/AKD16	1
Eimermacher / Sanacorp	Base cream DAC	86440517	1
Euro OTC	Base cream DAC	86490517	3
Eimermacher	Base cream DAC	86440517	1
Euro OTC	Base cream DAC	91310916	3
Euro OTC	Base cream DAC	91310616	1
EuRhoNoweda 11.3..17	Base cream DAC	91350916	2
Noweda	Base cream DAC	91340916	1
Euro OTC	Base cream DAC	91340916	2
Euro OTC	Base cream DAC	91330916	1
Euro OTC	Base cream DAC	91350916	2
Euro OTC	Base cream DAC	B55850317	1
Gehe	Base cream DAC	cdbgsdgb	1
Euro OTC	Base cream DAC	B72430816	1
Sanacorp	Base cream DAC		2
Bombastus	Base cream DAC	CHB289803	1
Fagron	Base cream DAC	T08-015887	1
Ensbona/ Sanacorp	Base cream DAC		2
/ Anzag	Base cream DAC		1
Bombastus	Base cream DAC	263423	1
Bombastus	Base cream DAC	265036	2
Bombastus	Base cream DAC	262001	1
Bombastus	Base cream DAC	263429	1
Bombastus	Base cream DAC	263291	1
Bombastus	Base cream DAC	265029	1
Bombastus	Base cream DAC	266637	3
Gehe	Base cream DAC	266525	1
Gehe	Base cream DAC	266663	1
Bombastus	Base cream DAC	266749	1
Bombastus	Base cream DAC	266752	1
Bombastus	Base cream DAC	266799	1
Bombastus	Base cream DAC	268815	2
Anzag 27.08.2012	Base cream DAC	268819	1
Bombastus	Base cream DAC	268821	1
Bombastus	Base cream DAC	268820	2
Bombastus	Base cream DAC	268925	6
PHOENIX am 08.06.2012	Base cream DAC	269997	1
Bombastus	Base cream DAC	269996	1
Bombastus	Base cream DAC	270000	3
Bombastus	Base cream DAC	272084	6
Bombastus	Base cream DAC	270157	1
Gehe	Base cream DAC	272000	1
Bombastus	Base cream DAC	271999	1
Bombastus	Base cream DAC	272001	1
Caelo	Base cream DAC	272084	2
Bombastus	Base cream DAC	272215	4
Bombastus	Base cream DAC	272091	1
Phoenix	Base cream DAC	272091	1
Bombastus	Base cream DAC	272086	2
Bombastus	Base cream DAC	272093	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream DAC	272219	1
Phoenix	Base cream DAC	272219	2
Phönix 22.11.12	Base cream DAC	272218	1
Bombastus	Base cream DAC	272216	1
Anzag 09.11.12	Base cream DAC	272216	1
Bombastus	Base cream DAC	272221	1
Bombastus	Base cream DAC	272464	3
Caelo	Base cream DAC	1081302	1
Caelo	Base cream DAC	1061303	1
Eimermacher	Base cream DAC	2770713	1
Caelo	Base cream DAC	213M-01309	1
Caelo	Base cream DAC	3121302	1
Hedinger	Base cream DAC	281718	1
Caelo	Base cream DAC	5091308	1
Gehe	Base cream DAC	4A131113	1
Caelo	Base cream DAC	513Q-01309	1
Euro OTC	Base cream DAC	80616-8	1
KEHR / Beyvers GmbH	Base cream DAC	71112PML01	1
Caelo	Base cream DAC	713I-01309	1
Caelo	Base cream DAC	8011302	1
Caelo	Base cream DAC	8111302	1
Caelo	Base cream DAC	913I-01309	1
Caelo	Base cream DAC	1004211309	1
Caelo	Base cream DAC	1013Q-01309	3
Noweda	Base cream DAC	11/30-T01	5
Euro OTC	Base cream DAC	11033-01	1
Caelo	Base cream DAC	11011405	1
EuRho	Base cream DAC	1108002-01	1
Euro OTC	Base cream DAC	1107040-01	1
Caelo	Base cream DAC	11087101	1
Euro OTC	Base cream DAC	1110037-01	2
Caelo	Base cream DAC	1109211309	1
Caelo	Base cream DAC	11101302	1
Euro OTC	Base cream DAC	11110023-01	1
Euro OTC	Base cream DAC	1111009-01	1
EuRho Anzag ek:33,16 03...	Base cream DAC	1111027-01	1
Europharm/Otto	Base cream DAC	1111028-01	2
Euro OTC	Base cream DAC	1111028-01	1
Phönix	Base cream DAC	1111029-01	1
Euro OTC	Base cream DAC	1111041-01	1
Euro OTC	Base cream DAC	1111042-01	1
EuRho	Base cream DAC	1111048-01	1
EuRho	Base cream DAC	1111050-01	1
Gehe	Base cream DAC	111105101	1
Euro OTC	Base cream DAC	1111051-01	2
Caelo	Base cream DAC	11206201	3
Caelo	Base cream DAC	11279001	2
Noweda	Base cream DAC	11229001	1
Caesar & Loretz GmbH, 4...	Base cream DAC	11278901	1
Caelo	Base cream DAC	11295701	1
Caelo	Base cream DAC	11278701	1
Caelo	Base cream DAC	11279401	1
Caelo	Base cream DAC	11295801	2
Caelo	Base cream DAC	11296601	1
Caelo	Base cream DAC	11296501	1
Caelo	Base cream DAC	11296301	1
Caelo	Base cream DAC	11296701	1
Caesar & Loretz	Base cream DAC	11296701	1
Caelo	Base cream DAC	11296901	1
Caelo	Base cream DAC	11297201	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	11297001	1
Fagron	Base cream DAC	11297301	1
Kehr	Base cream DAC	1130501	1
Caelo	Base cream DAC	11297301	1
Fagron	Base cream DAC	1130-T01	1
Caelo	Base cream DAC	11330401	1
Caelo	Base cream DAC	11330101	2
Caelo	Base cream DAC	11330201	1
Caelo	Base cream DAC	11330301	1
Caelo	Base cream DAC	11330601	1
Fagron	Base cream DAC	11B11-T05	4
Noweda	Base cream DAC	11K06-T01	5
Fagron	Base cream DAC	11K06-T01	1
Noweda	Base cream DAC	11130-T01	1
Noweda	Base cream DAC	11130-t01	1
Caelo	Base cream DAC	12000101	1
Caelo	Base cream DAC	12000401	2
Caelo	Base cream DAC	12000501	3
Caelo	Base cream DAC	12004901	2
Caelo	Base cream DAC	12005101	3
Caelo	Base cream DAC	12000601	1
Caelo	Base cream DAC	12005201	2
Caelo	Base cream DAC	12005501	5
Caelo	Base cream DAC	12005401	1
Noweda	Base cream DAC	12005601	1
Kehr	Base cream DAC	12005601	1
Caelo	Base cream DAC	12005701	6
Sanacorp	Base cream DAC	1202001	3
Doroderm GmbH/ Noweda	Base cream DAC	1202001	4
Caelo	Base cream DAC	1202001	2
doroderm GmbH	Base cream DAC	1202001	2
Euro OTC	Base cream DAC	1204004-01	5
Euro OTC	Base cream DAC	1204024-01	1
EuRho	Base cream DAC	1204022-01	1
Euro OTC	Base cream DAC	1204025-01	2
Euro OTC	Base cream DAC	1204023-01	1
Klenk	Base cream DAC	1205004-01	1
Caelo	Base cream DAC	1205007-01	2
Anzag 25.10.12 ek:33,15EUR	Base cream DAC	1205005-01	1
Euro OTC	Base cream DAC	1205018-01	3
Euro OTC	Base cream DAC	1205019-01	1
Euro OTC	Base cream DAC	1205040-01	1
Euro OTC	Base cream DAC	1205043-01	1
Euro OTC	Base cream DAC	1205044-01	2
Euro OTC	Base cream DAC	1205041-01	1
Caelo	Base cream DAC	12061701	1
Caelo	Base cream DAC	12061601	1
DAC 2011	Base cream DAC	12061901	1
Caelo	Base cream DAC	12062001	2
Caelo	Base cream DAC	12062601	2
Caesar & Loretz GmbH, 4...	Base cream DAC	12062601	1
Caelo	Base cream DAC	12062801	2
Caelo	Base cream DAC	12062701	2
Caelo	Base cream DAC	12063001	3
Euro OTC	Base cream DAC	1209037-01	4
Caelo	Base cream DAC	12062901	1
Caelo	Base cream DAC	12063101	1
Caelo	Base cream DAC	12071202	1
EuRho	Base cream DAC	1209037-01	1
Euro OTC	Base cream DAC	120903701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Base cream DAC	1209038-01	2
Anzag 14.01.2013	Base cream DAC	1210011-01	1
EuRho	Base cream DAC	1210011-01	2
Euro OTC	Base cream DAC	1210011-01	3
Bombastus	Base cream DAC	1210043-01	1
Euro OTC	Base cream DAC	1210044-01	4
Euro OTC	Base cream DAC	121005401	1
Euro OTC	Base cream DAC	1210053-01	1
Euro OTC	Base cream DAC	1210054-01	1
Hirsch	Base cream DAC	1210054-01	1
Euro OTC	Base cream DAC	1210055-01	3
Euro OCT	Base cream DAC	1210055-01	1
Euro OTC	Base cream DAC	1212025-01	1
Caelo	Base cream DAC	12101203	1
Caelo	Base cream DAC	12111202	1
Euro OTC	Base cream DAC	1212028-01	2
Caelo	Base cream DAC	12116001	1
anzag 29.01.13 ek: 33,15EUR	Base cream DAC	1212027-01	1
AHC 27.06.13 EK 32,95	Base cream DAC	1212025-01	1
Euro OTC	Base cream DAC	121202801	1
Euro OTC	Base cream DAC	121209-01	1
Euro OTC	Base cream DAC	1212029-01	3
Caelo	Base cream DAC	1213M-01309	1
Caelo	Base cream DAC	12159301	2
Caelo	Base cream DAC	12159201	1
Caelo	Base cream DAC	12159401	2
Caelo	Base cream DAC	12159601	5
Caelo	Base cream DAC	12167901	12
Caelo	Base cream DAC	12168001	2
Caelo	Base cream DAC	12168101	2
Kehr	Base cream DAC	12168301	1
Caelo	Base cream DAC	12168201	1
Caelo	Base cream DAC	12168301	1
Caelo	Base cream DAC	12174901	3
Caelo	Base cream DAC	12174801	3
Caelo	Base cream DAC	12175001	5
Finze	Base cream DAC	12174901	1
Caelo	Base cream DAC	12175201	2
Caelo/Gehe	Base cream DAC	12175101	1
Caelo	Base cream DAC	12175101	1
Caelo	Base cream DAC	12178301	4
Caelo	Base cream DAC	12178401	3
Caelo	Base cream DAC	12191701	3
Caelo	Base cream DAC	12191801	5
Sanacorp	Base cream DAC	12191801	1
Infecto/Phönix	Base cream DAC	12191901	1
Caelo	Base cream DAC	12192001	2
Caelo	Base cream DAC	12191902	1
Caelo	Base cream DAC	12192101	1
Caelo	Base cream DAC	12192201	2
Caelo	Base cream DAC	12240801	2
Noweda	Base cream DAC	12240801	1
Caelo	Base cream DAC	12240901	5
Hersteller: Beyvers	Base cream DAC	122/1	1
Caelo	Base cream DAC	12241001	5
Caelo	Base cream DAC	12241201	1
Caelo	Base cream DAC	12241101	1
phönix, 13.11.2012	Base cream DAC	12241201	1
Caelo	Base cream DAC	12245001	2
anzag, 20.11.2012	Base cream DAC	12245001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	12245101	2
Erika apotheke	Base cream DAC	1224501	1
Caelo	Base cream DAC	12249401	3
Caelo	Base cream DAC	12245201	1
Caelo	Base cream DAC	12249401	1
Caelo	Base cream DAC	12249701	2
Caelo	Base cream DAC	12262402	3
Sanacorp	Base cream DAC	12262402	1
Caelo	Base cream DAC	12262501	5
Caelo	Base cream DAC	12262601	3
Caelo	Base cream DAC	12262801	6
Caelo	Base cream DAC	1226290	1
Sanacorp	Base cream DAC	12262901	1
Caelo	Base cream DAC	12262901	3
Sanacorp	Base cream DAC	12262801	1
Phönix	Base cream DAC	12262901	1
Gehe	Base cream DAC	12262901	1
Sanacorp	Base cream DAC	12263001	1
Caelo	Base cream DAC	12263001	1
Kehr	Base cream DAC	12263101	1
Caelo	Base cream DAC	12263101	4
Caelo	Base cream DAC	12263201	6
Caelo	Base cream DAC	12263301	3
Phönix 01.02.13	Base cream DAC	12263201	1
Caelo	Base cream DAC	12263401	5
Caelo	Base cream DAC	12289701	4
Noweda	Base cream DAC	12263501	1
Sanacorp	Base cream DAC	12289701	1
Caelo	Base cream DAC	12303901	3
Sanacorp	Base cream DAC	12304701	1
Caelo	Base cream DAC	12304701	10
Caelo	Base cream DAC	12304801	9
Kehr	Base cream DAC	12304801	1
Sanacorp	Base cream DAC	12304801	1
Phönix	Base cream DAC	12304801	3
Caelo	Base cream DAC	12304901	3
Sanacorp	Base cream DAC	12304901	2
Caelo	Base cream DAC	12305001	5
Caelo	Base cream DAC	12305301	4
Caelo	Base cream DAC	12305101	1
Alliance Healthcare	Base cream DAC	12305301	1
Caelo	Base cream DAC	12339001	2
Caelo	Base cream DAC	12339101	3
Caesar & Loretz GmbH, 4...	Base cream DAC	12339101	1
Caelo	Base cream DAC	12339201	9
Noweda	Base cream DAC	12339201	1
Sanacorp	Base cream DAC	12339001	1
Krieger	Base cream DAC	12339301	1
Caelo	Base cream DAC	12339301	4
Euro OTC	Base cream DAC	12339301	1
Caelo	Base cream DAC	12339401	6
Caelo	Base cream DAC	12343801	6
Caelo	Base cream DAC	12347101	2
Gehe	Base cream DAC	12347301	1
Caelo	Base cream DAC	12347301	5
Caelo	Base cream DAC	12347401	5
Caelo	Base cream DAC	12384002	2
Caelo	Base cream DAC	1239201	1
Fagron	Base cream DAC	12C02-T08	5
Caelo	Base cream DAC	125044-1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Base cream DAC	129037-01	1
Fagron	Base cream DAC	12C26-T04	1
Fagron	Base cream DAC	12C02-T09	1
Fagron	Base cream DAC	12C02-T08	1
Großhandel	Base cream DAC	12C26-T08	1
Fagron	Base cream DAC	12C26-T08	2
Fagron	Base cream DAC	12G25-T03	3
Fagron	Base cream DAC	12G25T03	1
Fagron	Base cream DAC	12C26-T14	1
Fagron	Base cream DAC	12G25-T06	2
Fagron	Base cream DAC	12G25-T04	3
Noweda	Base cream DAC	12G25-T04	1
Caelo	Base cream DAC	130000601	1
Fagron	Base cream DAC	12G25-T10	1
Fagron	Base cream DAC	12G25-T12	2
Phönix	Base cream DAC	13000601	1
Caelo	Base cream DAC	13000701	2
Caelo	Base cream DAC	13000601	1
Caelo	Base cream DAC	13000501	1
Caelo	Base cream DAC	13000801	2
Phönix	Base cream DAC	13000701	1
Caelo	Base cream DAC	13007001	5
Sanacorp	Base cream DAC	13000901	1
Phönix	Base cream DAC	13007101	1
Caelo	Base cream DAC	13000901	1
Caelo	Base cream DAC	13007101	6
Caelo	Base cream DAC	13007201	3
phönix	Base cream DAC	13007101	1
Caelo	Base cream DAC	13007301	2
EuRho	Base cream DAC	1301022-01	2
Euro OTC	Base cream DAC	130102301	1
EuRho	Base cream DAC	1301023-01	1
Euro OTC	Base cream DAC	1301042-01	1
Euro OTC	Base cream DAC	1301040-01	1
Euro OTC	Base cream DAC	1301044-01	1
Caelo	Base cream DAC	13012401	10
Euro OTC	Base cream DAC	1301043-01	1
calo	Base cream DAC	13012401	1
Caelo	Base cream DAC	13012501	3
Caelo	Base cream DAC	13016801	2
Caelo	Base cream DAC	13012601	2
Caelo	Base cream DAC	1301681	1
Caelo	Base cream DAC	13016901	3
Caelo	Base cream DAC	13017001	3
Caelo	Base cream DAC	13017101	1
EuRho	Base cream DAC	1302009-01	1
Euro OTC	Base cream DAC	130201001	2
Euro OTC	Base cream DAC	1302036-01	1
EuRho	Base cream DAC	1302036-01	1
Euro OTC	Base cream DAC	1302037-01	2
Euro OTC	Base cream DAC	1304004-01	2
Caelo	Base cream DAC	13022401	3
EuRho	Base cream DAC	130400501	1
Euro OTC	Base cream DAC	1304005-01	2
Euro OTC	Base cream DAC	1304006-01	5
Euro OTC	Base cream DAC	1305012-01	2
Euro OTC	Base cream DAC	1305013-01	1
Euro OTC	Base cream DAC	1306001-01	6
Caelo	Base cream DAC	1306002-01	1
Euro OTC	Base cream DAC	1306002-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro/Noweda	Base cream DAC	1306008-01	3
Euro OTC	Base cream DAC	1306007-01	1
Euro OTC	Base cream DAC	1306009-01	2
Euro OTC	Base cream DAC	1306020-01	2
Euro OTC	Base cream DAC	1306031-01	1
Euro OTC	Base cream DAC	1306019-01	2
Caelo	Base cream DAC	1306031-01	1
Euro OTC	Base cream DAC	130603001	1
Euro OTC	Base cream DAC	1306042-01	1
Euro OTC	Base cream DAC	1306044-01	2
Gehe	Base cream DAC	1306044-01	1
Caelo	Base cream DAC	1307031-01	4
Euro OTC	Base cream DAC	130703101	1
Euro OTC	Base cream DAC	130703601	2
Caelo	Base cream DAC	130704501	1
Euro OTC	Base cream DAC	1307046-01	1
Euro OTC	Base cream DAC	1307035-01	1
Euro OTC	Base cream DAC	130704801	1
Euro OTC	Base cream DAC	1308011-01	1
Euro OTC	Base cream DAC	130802-01	1
Euro OCT	Base cream DAC	1307047-01	1
Caelo	Base cream DAC	13082101	10
Caelo	Base cream DAC	13082203	9
Caelo	Base cream DAC	13083101	3
Euro OTC	Base cream DAC	1309005-01	3
Euro OTC	Base cream DAC	1311007-01	3
Euro OTC	Base cream DAC	131100901	1
Euro OTC	Base cream DAC	1311009-01	1
Euro OTC	Base cream DAC	1311018-01	1
Caelo	Base cream DAC	13111101	3
Euro OTC	Base cream DAC	1311031-01	1
Caelo	Base cream DAC	13111201	4
Caelo	Base cream DAC	13111301	5
Caelo	Base cream DAC	13111401	2
Phönix	Base cream DAC	13111401	2
Euro OTC	Base cream DAC	13113101	1
Caelo	Base cream DAC	13115601	8
Caelo	Base cream DAC	13115701	3
Hedinger/Phönix	Base cream DAC	13115801	3
Caelo	Base cream DAC	13115801	2
Caelo	Base cream DAC	13115901	3
Caelo	Base cream DAC	13116001	2
Caelo	Base cream DAC	13119901	1
Caelo	Base cream DAC	13120001	2
Euro OTC	Base cream DAC	1312001-01	2
Caelo	Base cream DAC	13120101	5
Caelo	Base cream DAC	1312001-01	1
Caelo	Base cream DAC	13120201	4
Caelo	Base cream DAC	13120301	7
Caelo	Base cream DAC	13125001	4
Caelo	Base cream DAC	13125101	2
Caelo	Base cream DAC	13139702	10
Caesar & Loretz GmbH, 4...	Base cream DAC	13139702	1
Noweda	Base cream DAC	13139702	1
Alliance Healthcare	Base cream DAC	13139801	1
Noweda	Base cream DAC	13139901	1
Caesar /phönix	Base cream DAC	13139901	2
Caelo	Base cream DAC	13140001	5
Caelo	Base cream DAC	13158401	1
Caelo	Base cream DAC	1315601	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	13158501	9
Caelo	Base cream DAC	13158440	1
Caelo	Base cream DAC	13158601	4
Caelo	Base cream DAC	13158801	3
Caelo	Base cream DAC	13173301	2
Caelo	Base cream DAC	13173601	14
Herstellung	Base cream DAC	13173601	1
Sanacorp	Base cream DAC	13173601	1
Caelo	Base cream DAC	13173701	5
Caelo	Base cream DAC	1318701	1
Sanacorp	Base cream DAC	13173701	1
Caelo	Base cream DAC	13187801	15
Caelo	Base cream DAC	13187701	2
Caelo	Base cream DAC	13187901	10
Caelo	Base cream DAC	13188001	7
Alliance Healthcare	Base cream DAC	13188101	1
Caelo	Base cream DAC	13188101	3
Caelo	Base cream DAC	13188401	7
DAC 2011	Base cream DAC	13188401	1
Caelo	Base cream DAC	13188501	9
Caelo	Base cream DAC	13188601	4
Caelo	Base cream DAC	13188701	11
Sanacorp	Base cream DAC	13188801	1
Caelo	Base cream DAC	13188801	3
Caelo	Base cream DAC	1319702	1
Caelo	Base cream DAC	13198401	1
Caelo	Base cream DAC	13198501	3
Gehe	Base cream DAC	13198501	1
Caelo	Base cream DAC	13198601	8
Noweda	Base cream DAC	13198601	2
Caelo	Base cream DAC	13198801	20
Caelo	Base cream DAC	13198901	8
Caelo	Base cream DAC	13199001	11
Caelo	Base cream DAC	13199101	19
Sanacorp	Base cream DAC	13199101	3
Euro OTC	Base cream DAC	13201101	1
Caelo	Base cream DAC	132011-01	1
Gehe	Base cream DAC	13199101	1
Caelo	Base cream DAC	13249702	3
Caelo	Base cream DAC	13249801	6
Fagron	Base cream DAC	13249901	1
Caelo	Base cream DAC	13249901	2
Sanacorp	Base cream DAC	13250001	1
Caelo	Base cream DAC	13250001	8
Caelo	Base cream DAC	13252101	5
Caelo	Base cream DAC	13310901	7
Gehe	Base cream DAC	13310901	1
Caelo	Base cream DAC	13311201	10
Caelo	Base cream DAC	13311101	3
Caelo	Base cream DAC	13311301	11
Sanacorp	Base cream DAC	13311301	2
Caelo	Base cream DAC	13311601	4
Caelo	Base cream DAC	13311502	1
Caelo	Base cream DAC	13311606	1
Caelo	Base cream DAC	13311701	5
Fagron	Base cream DAC	13311701	1
Phönix	Base cream DAC	13311801	1
Sanacorp	Base cream DAC	13312001	1
Caelo	Base cream DAC	1331201	1
Caelo	Base cream DAC	13312101	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream DAC	13312201	5
Caelo	Base cream DAC	13312301	4
Sanacorp	Base cream DAC	13312301	1
Caelo	Base cream DAC	13312401	2
Caelo	Base cream DAC	1331801	1
Sanacorp	Base cream DAC	13377201	1
celo/phönix	Base cream DAC	13377901	3
Fagron	Base cream DAC	13378001	2
Sanacorp	Base cream DAC	13380901	1
Sanacorp	Base cream DAC	13381001	1
Sanacorp	Base cream DAC	13381101	1
D'ORODERM	Base cream DAC	1346001	1
doroderm GmbH	Base cream DAC	1346001	4
Caelo	Base cream DAC	1387801	1
Caelo	Base cream DAC	13812013	1
Fagron	Base cream DAC	13A08-T06-000034	1
Fagron	Base cream DAC	13A09T10000232	1
Fagron	Base cream DAC	13A09-T10-000232	3
Caelo	Base cream DAC	13A09-T10-232	2
Fagron	Base cream DAC	13A09-T12-000238	1
Fagron	Base cream DAC	13A09-T12-236	1
Fagron	Base cream DAC	13A09-T14-000244	1
Fagron	Base cream DAC	13A09-T14-000244	1
Fagron	Base cream DAC	13A09-T18-000259	2
Fagron	Base cream DAC	13A09-T19-00194	1
Fagron	Base cream DAC	13C05-T-003297	1
Fagron	Base cream DAC	13C05-T03-001755	1
Fagron	Base cream DAC	13C05-T03-00175	1
Eurho/Phoenix	Base cream DAC	13C05-T03-001756	1
Fagron	Base cream DAC	13C05-T03-001756	8
Fagron	Base cream DAC	13C05-T04-001760	1
Fagron	Base cream DAC	13C05-T06-001767	1
Fagron	Base cream DAC	13C05-T09-001778	3
Fagron	Base cream DAC	13C05-T10-003297	1
Fagron	Base cream DAC	13C05-T10-003301	1
Fagron	Base cream DAC	13E21-T04-003893	1
Fagron	Base cream DAC	13E20-T01-003865	1
Fagron	Base cream DAC	13E21-T05-003897	7
Phönix	Base cream DAC	13E21-T05-003897	1
Caelo	Base cream DAC	13E21-T06-003900	1
Fagron	Base cream DAC	13E21-T06-003900	1
Caelo	Base cream DAC	13E21-T08-003906	2
Fagron	Base cream DAC	13E22-T04-003975	1
Fagron	Base cream DAC	13E22-T05-003978	1
Fagron	Base cream DAC	13E22-T07-003984	3
Großhandel	Base cream DAC	13E21-T09-003911	1
Gehe	Base cream DAC	13E22T05003978	1
Fagron	Base cream DAC	13E22T09003991	1
Fagron	Base cream DAC	13E23-004023	1
Fagron	Base cream DAC	13E22-T08-003987	2
Fagron	Base cream DAC	13E23-T06-004008	3
Fagron	Base cream DAC	13E23-T07-004011	4
Fagron	Base cream DAC	13E23-T08-004015	1
Fagron	Base cream DAC	13E23-T09	3
Fagron	Base cream DAC	13h14t02006207	1
Caelo	Base cream DAC	13E23-T10-004023	2
Fagron	Base cream DAC	13E23-T10-004023	1
Fagron	Base cream DAC	13E23-T12	1
Caelo	Base cream DAC	13H14-T01-006203	2
Fagron	Base cream DAC	13H14T02006207	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Base cream DAC	13H14-T02-006207	1
Fagron	Base cream DAC	13H14-T01-006202	1
Caelo	Base cream DAC	13H14T05-006551	1
Maros/Phönix	Base cream DAC	13H14-To3-006212	1
Phoenix	Base cream DAC	14036001	1
Caelo	Base cream DAC	1409101	1
Caelo	Base cream DAC	14090501	1
Caelo	Base cream DAC	14121-01309	1
Caelo	Base cream DAC	14144202	1
Caelo	Base cream DAC	144144201	1
p.w beyers gmbh,phönix	Base cream DAC		4
P.W.Beyvers GmbH	Base cream DAC		3
Hersteller: Beyvers	Base cream DAC		1
Caelo	Base cream DAC	14A10-T06-009259	1
Fagron	Base cream DAC	14D08-T03-01168	1
Euro OTC	Base cream DAC	150200801	1
Caelo	Base cream DAC	14D08T04-011692	1
Phönix	Base cream DAC	15272601	1
Euro OTC	Base cream DAC	1520202801	1
Phönix	Base cream DAC	15450601	1
Caelo	Base cream DAC	160704	1
Celo	Base cream DAC	16005101	1
Caelo	Base cream DAC	16B290713	1
Caelo	Base cream DAC	17061305	1
Fagron	Base cream DAC	16C14T03-031205	1
Caelo	Base cream DAC	1713M-01309	1
Caelo	Base cream DAC	19101302	2
Caelo	Base cream DAC	18071302	1
P.W.Beyvers GmbH	Base cream DAC	207/1	2
Caelo	Base cream DAC	19021304	1
Gehe	Base cream DAC	207/1	4
Caelo	Base cream DAC	2025-T05	1
Gehe	Base cream DAC	207/2	1
Anzag	Base cream DAC	207/2	1
Ferdi	Base cream DAC	21220113	1
P.W.Beyvers GmbH,Phönix	Base cream DAC	207/2	1
PHOENIX am 04.06.2012	Base cream DAC	219/1	3
P.W.Beyvers GmbH	Base cream DAC	219/1	1
P.W. Beyvers	Base cream DAC	219/1	2
Gehe	Base cream DAC	219/1	3
P.W.Beyvers GmbH	Base cream DAC	219/2	1
Gehe	Base cream DAC	219/2	3
Caelo	Base cream DAC	22031303	1
Caelo	Base cream DAC	219/2	2
Caelo	Base cream DAC	2213V01309	1
Caelo	Base cream DAC	22041306	1
P.W.Beyvers GmbH	Base cream DAC	229/1	1
P.W. Beyvers	Base cream DAC	229/1	1
P.W.Beyvers	Base cream DAC	229/1	2
P.W.Beyvers	Base cream DAC	229/1	1
Gehe	Base cream DAC	229/2	4
P.W.Beyvers GmbH	Base cream DAC	229/2	3
anzag, 05.11.2012	Base cream DAC	229/2	1
P.W.Beyvers GmbH / Phoe...	Base cream DAC	234/1	1
P.W.Beyvers	Base cream DAC	229/2	1
Beyers	Base cream DAC	234/1	1
Gehe	Base cream DAC	234/1	1
Gehe	Base cream DAC	234/2	4
Beyvers	Base cream DAC	234/2	1
P.W.Beyvers	Base cream DAC	234/1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
P.W. Beyvers	Base cream DAC	234/1	1
Gehe 17.01.13	Base cream DAC	234/2	1
P.W. Beyvers	Base cream DAC	234/3	2
P.W. Beyvers	Base cream DAC	234/2	1
P.W.Beyvers Gmbh,Phönix	Base cream DAC	234/3	2
Beyvers	Base cream DAC	234/3	1
P.W.Bayvers GmbH	Base cream DAC	234-2	1
Caelo	Base cream DAC	24091301	1
Ensbona/Eimermacher	Base cream DAC	24550213	1
P.W.Beyvers GmbH / Phoe...	Base cream DAC	243/3	1
Caelo	Base cream DAC	25011401	1
Caelo	Base cream DAC	26061204	1
Caelo	Base cream DAC	25021402	1
Gehe	Base cream DAC	272367	1
Caelo	Base cream DAC	272369	2
Gehe	Base cream DAC	272371	1
Gehe	Base cream DAC	272373	1
Bombastus	Base cream DAC	272375	2
Bombastus	Base cream DAC	275172	1
Bombastus	Base cream DAC	275173	1
Bombastus	Base cream DAC	272373	1
Gehe	Base cream DAC	275147	1
Bombastus	Base cream DAC	275181	2
Bombastus	Base cream DAC	275182	1
Caelo	Base cream DAC	275184	2
Bombastus	Base cream DAC	275277	1
Bombastus	Base cream DAC	275273	1
Bombastus	Base cream DAC	275279	1
Bombastus	Base cream DAC	275280	1
Bombastus	Base cream DAC	275438	6
Bombastus	Base cream DAC	275282	1
Caelo	Base cream DAC	275438	2
AHD, 05.04.2013	Base cream DAC	275439	1
Bombastus	Base cream DAC	275440	1
Bombastus	Base cream DAC	275590	2
Bombastus	Base cream DAC	275443	1
Anzag 19.02.2013	Base cream DAC	275444	1
Bombastus	Base cream DAC	275447	1
Bombastus	Base cream DAC	275445	1
Caelo	Base cream DAC	275591	1
Bombastus	Base cream DAC	275591	1
Phoenix	Base cream DAC	275591	1
Bombastus	Base cream DAC	276944	2
Bombastus	Base cream DAC	276948	1
Bombastus	Base cream DAC	276950	2
Bombastus	Base cream DAC	277099	14
phönix	Base cream DAC	277099	1
Bombastus	Base cream DAC	277112	1
Gehe	Base cream DAC	277109	1
Caelo	Base cream DAC	277111	1
Bombastus	Base cream DAC	277248	1
Bombastus	Base cream DAC	277110	1
Gehe	Base cream DAC	278392	2
Bombastus	Base cream DAC	278396	2
Bombastus	Base cream DAC	278397	1
Bombastus	Base cream DAC	278470	2
Gehe	Base cream DAC	278401	1
Apomix / 11	Base cream DAC	278472	2
Bombastus	Base cream DAC	278472	2
Bombastus	Base cream DAC	278474	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream DAC	278475	1
Bombastus	Base cream DAC	278476	1
Bombastus	Base cream DAC	278478	1
Bombastus	Base cream DAC	278479	2
Phönix	Base cream DAC	278475	1
Bombastus	Base cream DAC	278553	1
Bombastus	Base cream DAC	278555	4
Bombastus	Base cream DAC	278556	7
Bombastus	Base cream DAC	278557	2
Bombastus	Base cream DAC	278560	2
Bombastus	Base cream DAC	278558	1
Bombastus	Base cream DAC	279829	1
Bombastus	Base cream DAC	278562	2
Bombastus	Base cream DAC	278561	1
Bombastus	Base cream DAC	278559	1
Bombastus	Base cream DAC	279958	3
Caelo	Base cream DAC	279963	1
Bombastus	Base cream DAC	279960	1
Phönix, 28.10.2013	Base cream DAC	279961	1
Bombastus	Base cream DAC	279964	3
Bombastus	Base cream DAC	279961	1
Bombastus	Base cream DAC	280077	1
Holdermann	Base cream DAC	280078	1
AHD, 02.12.13	Base cream DAC	280079	1
Bombastus	Base cream DAC	280079	2
Bombastus	Base cream DAC	280080	4
Gehe	Base cream DAC	280082	1
Caelo	Base cream DAC	28031401	1
Caelo	Base cream DAC	28031403	1
Caelo	Base cream DAC	280913B	1
Caelo	Base cream DAC	28031404	1
Bombastus	Base cream DAC	281239	1
gehe	Base cream DAC	281235	1
Gehe	Base cream DAC	281237	1
Bombastus	Base cream DAC	281240	1
Apomix / AHD	Base cream DAC	281328006/14	1
Bombastus	Base cream DAC	281328	3
Bombastus	Base cream DAC	281329	3
Bombastus	Base cream DAC	281331	4
Bombastus	Base cream DAC	281332	1
Bombastus	Base cream DAC	281486	3
Bombastus	Base cream DAC	281335	1
Bombastus	Base cream DAC	282723	2
Gehe	Base cream DAC	282769	1
Sanacorp	Base cream DAC	282772	1
Bombastus	Base cream DAC	282771	1
Eimermacher	Base cream DAC	28310913	2
Bombastus	Base cream DAC	282778	1
Bombastus	Base cream DAC	283917	1
Gehe	Base cream DAC	283975	1
Sanacorp	Base cream DAC	283972	2
Bombastus	Base cream DAC	288058	1
Bombastus	Base cream DAC	286086	2
Caelo	Base cream DAC	2901311309	1
Caelo	Base cream DAC	30111201	1
Bombastus	Base cream DAC	299325	1
Caelo	Base cream DAC	3013M-01309	1
P.W.Beyvers	Base cream DAC		1
Noweda	Base cream DAC	30241211	4
Caelo	Base cream DAC	31081205	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
P.W. Beyvers	Base cream DAC		1
Beyvers	Base cream DAC		3
Beyvers	Base cream DAC	319/1	2
Gehe	Base cream DAC	319/1	2
P.W.Beyers/Phoenix	Base cream DAC	319/1	1
gehe 24.07.2013	Base cream DAC	319/1	1
P. W. Beyvers	Base cream DAC	319/1	1
P.W.Beyers/Phoenix	Base cream DAC	319/2	1
P.W. Beyvers	Base cream DAC	329/1	2
Beyvers	Base cream DAC	319/2	1
Gehe	Base cream DAC	329/1	3
Beyvers	Base cream DAC	319-1	1
Beyvers	Base cream DAC	329/1	1
Beyvers GmbH/ Noweda	Base cream DAC	329/1	1
P.W.Beyvers GmbH,Phönix	Base cream DAC	329/1	1
P.W. Beyvers	Base cream DAC	329/2	4
Gehe	Base cream DAC	329/2	1
Beyvers GmbH	Base cream DAC	329/2	1
P. W. Beyvers GmbH	Base cream DAC	329/2	3
Beyvers	Base cream DAC	329/2	2
Gehe	Base cream DAC	340/1	2
Beyvers/Anzag	Base cream DAC	340/1	1
P.W.Weyers GmbH	Base cream DAC	340/1	1
gehe 28.11.13	Base cream DAC	340/2	1
P.W.Beyvers	Base cream DAC	340/1	1
Beyvers	Base cream DAC	340/2	3
Beyvers	Base cream DAC	340/1	1
P.W. Beyvers / Anzag	Base cream DAC	340/2	1
P.W.Beyvers	Base cream DAC	340/2	3
P.W.Beyvers GmbH/Anzag	Base cream DAC	340/2	1
P.W.Beyvers/Jenne 18.12...	Base cream DAC	340/2	1
Beyvers/Gumo	Base cream DAC	340/2	1
P.W.Beyvers GmbH/Phönix	Base cream DAC	340/2	1
Caelo	Base cream DAC	3413E-01309	1
Caelo	Base cream DAC	3513I-01309	1
PW Beyves Alliance	Base cream DAC	412/3	1
Fagron	Base cream DAC	3613E-01309	1
Eimermacher	Base cream DAC	36140913	1
Caelo	Base cream DAC	4113E-01309	1
Caelo	Base cream DAC	4121	1
anzag, 20.11.2012	Base cream DAC	412245001	1
P.W.Beyvers GmbH,Phönix	Base cream DAC	440/1	2
P.W.Beyvers GmbH,Phönix	Base cream DAC	440/2	1
Caelo	Base cream DAC	4913E01309	1
Caelo	Base cream DAC	4512Q-01309	1
Noweda	Base cream DAC	5099500	3
Caelo	Base cream DAC	5012Q-01309	1
Noweda	Base cream DAC	50995000	1
Henry Lamotte	Base cream DAC	5100000	1
Lamotte	Base cream DAC	5099700	1
Noweda	Base cream DAC	5099900	1
Lamotte	Base cream DAC	5099800	1
Henry Lamotte Oils GmbH	Base cream DAC	5151700	1
Henry Lamotte	Base cream DAC	5163400	1
Noweda	Base cream DAC	5100000	1
Phönix Gotha	Base cream DAC	5163300	1
Henry Lamotte Oils GmbH	Base cream DAC	5163300	1
Lamotte	Base cream DAC	5163400	1
Henry Lamotte Oils GmbH	Base cream DAC	5215300	3
Noweda	Base cream DAC	5163300	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte	Base cream DAC	5215400	1
Phönix Gotha	Base cream DAC	5215600	2
Henry Lamotte Oils GmbH	Base cream DAC	5215500	3
Henry Lamotte Oils	Base cream DAC	5215600	1
Henry Lamotte Oils GmbH	Base cream DAC	5215700	1
Phönix Gotha	Base cream DAC	5215800	1
Lamotte/Otto	Base cream DAC	5215800	2
Noweda	Base cream DAC	52701012	1
Caelo	Base cream DAC	5289100	1
Ensbona/Eimermacher	Base cream DAC	52701012	1
Phoenix	Base cream DAC	5289100	1
Henry Lamotte	Base cream DAC	5289000	1
Lamotte/Otto	Base cream DAC	5321000	1
Anzag/Henry Lamotte	Base cream DAC	5289100	1
Henry Lamotte/Anzag	Base cream DAC	5289200	1
Henry Lamotte Oils GmbH	Base cream DAC	5321000	5
Phönix Gotha	Base cream DAC	5321000	1
Lamotte 001	Base cream DAC	5321000	1
Henry Lamotte	Base cream DAC	5321000	1
Henry Lamotte Oils/ Gehe	Base cream DAC	5321000	2
Anzag 11.03.13	Base cream DAC	5321100	1
Henry Lamotte	Base cream DAC	5321100	1
Henry Lamotte Oils	Base cream DAC	5321200	1
Noweda/Lamotte Oils	Base cream DAC	5321200	1
Caelo	Base cream DAC	5321200	2
Phönix Gotha	Base cream DAC	5322100	1
Anzag 11.03.2013	Base cream DAC	5322200	1
Henry Lamotte Oils GmbH	Base cream DAC	5322100	1
Lamotte/Phoenix	Base cream DAC	5322200	1
Henry Lamotte/ Noweda	Base cream DAC	532220	1
Henry Lamotte Oils GmbH	Base cream DAC	5322200	1
Henry Lamotte Oils	Base cream DAC	5322400	1
Henry Lamotte	Base cream DAC	5322300	3
Lamotte 001	Base cream DAC	5322400	1
Alliance Healthcare	Base cream DAC	5322400	1
Caelo	Base cream DAC	538/2	2
Henry Lamotte Oils GmbH	Base cream DAC	5407700	5
post apotheke	Base cream DAC	540770	1
Henry Lamotte	Base cream DAC	5407700	1
lamotte	Base cream DAC	5407800	1
Lamotte/ Noweda	Base cream DAC	5407900	1
Lamotte/AHD	Base cream DAC	5407900	1
Lamotte / Anzag	Base cream DAC	5408000	1
Henry Lamotte	Base cream DAC	5408000	1
Caelo	Base cream DAC	5408000	1
Henry Lamotte	Base cream DAC	5408100	1
henry Lamotte / Noweda	Base cream DAC	5408000	1
Lamotte	Base cream DAC	5408000	1
Lamotte/Anzag	Base cream DAC	5408100	1
Lamotte / Anzag	Base cream DAC	5408200	1
Lamotte / Anzag	Base cream DAC	5408100	1
Henry Lamotte/Anzag	Base cream DAC	5408300	2
Anzag/Lamotte	Base cream DAC	5408300	1
phönix	Base cream DAC	5408200	1
Henry Lamotte	Base cream DAC	5408300	1
Henry Lamotte	Base cream DAC	5489800	1
Henry Lamotte Oils GmbH	Base cream DAC	5489800	1
sanacorp	Base cream DAC	548900	1
Henry Lamotte Oils GmbH...	Base cream DAC	5489900	1
Caelo	Base cream DAC	5490000	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils GmbH	Base cream DAC	5490100	1
Caelo	Base cream DAC	5490100	1
Alliance Health 17.01.2...	Base cream DAC	5490100	1
Henry Lamotte	Base cream DAC	5490100	1
Lamotte	Base cream DAC	5537000	1
henry lamotte oils gmbh	Base cream DAC	5537000	1
Lamotte/ Noweda	Base cream DAC	5537300	1
Lamotte	Base cream DAC	5650700	3
Caelo	Base cream DAC	5537400	1
Euro OTC	Base cream DAC	5650900	1
Lamotte	Base cream DAC	5785800	1
Lamotte	Base cream DAC	5787800	1
Lamotte	Base cream DAC	5787900	2
Lamotte	Base cream DAC	5651000	1
Lamotte	Base cream DAC	5788000	1
Lamotte	Base cream DAC	5788200	2
Lamotte	Base cream DAC	5909300	3
Henry Lamotte	Base cream DAC	5909300	1
Lamotte/Noweda	Base cream DAC	5909600	1
Lamotte	Base cream DAC	5909600	2
Lamotte/Noweda	Base cream DAC	6002237001	1
Lamotte	Base cream DAC	6002237002	1
henry lamotte	Base cream DAC	6002237003	1
Lamotte	Base cream DAC	59096006002237002	1
Henry Lamotte	Base cream DAC	6002237003	1
Caelo	Base cream DAC	6002237003	2
Henry Lamotte Oils GmbH	Base cream DAC	6002237004	1
Bombastus	Base cream DAC	6862441	1
Lamotte	Base cream DAC	6002237005	1
Lamotte	Base cream DAC	60260	1
ensbona	Base cream DAC	74000513	1
Eimermacher	Base cream DAC	97870214	2
Caelo	Base cream DAC	AKD6	1
Ensbona	Base cream DAC	98480713	1
Euro OTC	Base cream DAC	Ch1204005-01	1
Caelo	Base cream DAC	CH5490000	1
VDL;4,87EUR;12.06.15	Base cream DAC	14232301	2
Euro OTG/ Euro OTG Pharma	Base cream DAC	1501041-01	1
P&M Cosmetics GmbH & Co...	Base cream DAC	150317	1
Euro OTC	Base cream DAC	1604023-01	1
Caelo	Base cream DAC	160976012	1
Caelo	Base cream DAC	1623601	1
Caelo	Base cream DAC	16232801	1
Sanacorp, 8,04EUR, 24.1.17	Base cream DAC	16350101	1
Bombastus	Base cream DAC	305729	1
henry lamotte	Base cream DAC	5537100	1
henry lamotte	Base cream DAC	5650700	1
Henry Lamotte oils GmbH	Base cream DAC	5787800	2
Eimermacher Direktbezug	Base cream DAC	57080315	1
Henry Lamotte	Base cream DAC	5788200	1
Lamotte/AHD	Base cream DAC	6002237009	1
Lamotte/AHD	Base cream DAC	5909900	1
Henry Lamotte /AHD	Base cream DAC	5910000	1
Henry Lamotte /AHD	Base cream DAC	6003682017	1
Lamotte/AHD	Base cream DAC	6004860003	1
Euro OTC	Base cream DAC	6004860001	1
Lamotte/AHD	Base cream DAC	6003682005	1
Lamotte/AHD	Base cream DAC	8003325008	1
Bombastus	Base cream DAC	266524	1
Bombastus	Base cream DAC	269993	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Dorfmeister	Base cream DAC	1205041-01	1
Caelo	Base cream DAC	11279301	3
Anzag	Base cream DAC	1205006-01	1
Euro OTC	Base cream DAC	1205004-01	1
Euro OTC	Base cream DAC	1209039-01	1
Gehe	Base cream DAC	12005301	1
Euro OTC	Base cream DAC	1212027-01	1
GEHE	Base cream DAC	12174801	1
Caelo	Base cream DAC	1230901	1
Caelo	Base cream DAC	12347201	1
Bombastus	Base cream DAC	272372	1
Gehe	Base cream DAC	275277	1
Gehe	Base cream DAC	275590	1
Gehe	Base cream DAC	279961	1
Bombastus	Base cream DAC	275175	1
Bombastus	Base cream DAC	280082	2
Anzag	Base cream DAC	5215600	1
Gehe	Base cream DAC	287941	1
Bombastus	Base cream DAC	282774	1
henry lamotte	Base cream DAC	5489800	2
Fagron	Base cream DAC	5321100	1
Anzag	Base cream DAC	5321000	1
henry lamotte	Base cream DAC	5490100	1
Caelo	Base cream DAC	1306020-01	5
Caelo	Base cream DAC	131203801	1
Gutenberg-Apotheke	Base cream DAC	13378201	1
DORODERM GmbH & Co.KG	Base cream DAC	1346001	2
Caelo	Base cream DAC	14053201	1
Gutenberg-Apotheke	Base cream DAC	14035801	1
Fagron	Base cream DAC	14223201	1
Fagron	Base cream DAC	14F05-T07-012966	1
Bombastus	Base cream DAC	287193	1
Bombastus	Base cream DAC	289591	1
Bombastus	Base cream DAC	292258	1
Bombastus	Base cream DAC	294837	1
Henry Lamotte oils GmbH	Base cream DAC	5490100	1
Caelo	Base cream DAC	5537300	1
Henry Lamotte oils GmbH	Base cream DAC	5650900	1
Henry Lamotte oils GmbH	Base cream DAC	5650600	1
Henry Lamotte/Anzag	Base cream DAC	5651100	1
Lamotte/AHD	Base cream DAC	5651300	1
Henry Lamotte oils GmbH	Base cream DAC	5651100	1
Anzah / Henry lamotte	Base cream DAC	57788600	1
Caelo	Base cream DAC	12344001	1
Caelo	Base cream DAC	1301023-01	1
Sanacorp	Base cream DAC	13007201	1
Sanacorp	Base cream DAC	13007301	1
Caelo	Base cream DAC	13012301	2
Sanacorp	Base cream DAC	13082101	1
Sanacorp	Base cream DAC	13082203	1
Fagron	Base cream DAC	13111101	1
Sanacorp	Base cream DAC	13158401	2
Sanacorp	Base cream DAC	13139801	1
Henry Lamotte/Sanacorp	Base cream DAC	13140001	1
Sanacorp	Base cream DAC	13187901	1
Sanacorp	Base cream DAC	13188701	1
Sanacorp	Base cream DAC	13310901	1
Sanacorp	Base cream DAC	1317301	1
Caelo	Base cream DAC	13311301-3462312170	1
Sanacorp	Base cream DAC	13347201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Base cream DAC	13312201	1
Fargon	Base cream DAC	13E20-T01-003865	1
Fargon	Base cream DAC	13A09-T19-001952	1
Fargon	Base cream DAC	13H14-T10-007314	1
Ichtyol	Base cream DAC	5407700	1
Lamotte	Base cream DAC	5322200	1
Bombastus	Base cream DAC	282887	1
Lamotte	Base cream DAC	5489900	1
Fargon	Base cream DAC	13A09-T19-001949	1

- 85 072 spectra from 500 *Apo-Ident* customers from a total of 24 898 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Base cream DAC* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Base cream DAC* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	620	0	34 924
Type B	0	5367	45	88 382
Type C	0	5845	435	85 065

The substance/substance group *Base cream DAC* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.0323 %)
Type B	100.0000 % (> 99.9763 %)	99.1685 % (> 99.1131 %)
Type C	100.0000 % (> 99.6638 %)	93.0732 % (> 93.0255 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30923	31252	9.04	14.69
31043	34030	5.70	21.47
31252	31252	0.00	15.63
31474	31474	0.00	16.37
31709	31709	0.00	20.40
31710	31710	0.00	20.07
31711	31711	0.00	21.27
31712	31712	0.00	22.51
31713	31713	0.00	21.52
31714	31714	0.00	21.48
31715	31715	0.00	22.54
31716	31716	0.00	21.06
31964	31964	0.00	21.37
32345	31716	4.27	14.92
34030	34030	0.00	20.44

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Basis Cordes® RK
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30892-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Basis Cordes® RK

Special notes

When selecting the *Basis Cordes® RK* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Basis Cordes® RK*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Ichthyol	Basis Cordes® RK	12A114	30892	60	1405097
Ichthyol	Basis Cordes® RK	13A038	31078	40	20130301
Ichthyol	Basis Cordes® RK	14A087	31660	40	1511145

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 140 spectra of 3 reference samples from the substance/substance group *Basis Cordes*[®] RK. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 404 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 130 spectra from 4 reference samples of the substance/substance group *Basis Cordes*[®] RK.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Ichthyol	Basis Cordes [®] RK	10A001	30225	30
Ichthyol	Basis Cordes [®] RK	14A087	31660 [†]	20
Ichthyol	Basis Cordes [®] RK	16A035	32951	40
Ichthyol	Basis Cordes [®] RK	16A034	33617	40

- 93 744 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 499 spectra from 194 *Apo-Ident* customers from 55 batches from the substance/substance group *Basis Cordes*[®] RK.
- These include spectra of independent samples from 52 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Ichthyol	Basis Cordes [®] RK	10071401	1
Ichthyol/Sanacorp	Basis Cordes [®] RK	12A167	1
Ichthyol / Noweda	Basis Cordes [®] RK	13A039	3
Ichthyol-Gesellschaft/P...	Basis Cordes [®] RK	13A039	1
Ichthyol	Basis Cordes [®] RK	13A039	2
Ichthyol/Noweda	Basis Cordes [®] RK	13A103	2
Ichthyol/Otto	Basis Cordes [®] RK	13A103	1
Ichthyol-Gesellschaft/N...	Basis Cordes [®] RK	13A103	1
Ichthyol-Gesellschaft	Basis Cordes [®] RK	13A039	2
Ichthyol-Gesellschaft/G...	Basis Cordes [®] RK	13A103	2
Caelo	Basis Cordes [®] RK	13A103	1
Ichthyol /Noweda	Basis Cordes [®] RK	13A103	3
Phönix Gotha	Basis Cordes [®] RK	13A103	2
Ichthyol/Gehe	Basis Cordes [®] RK	13A103	1
Ichtyol	Basis Cordes [®] RK	13A103	2
Ichthyol	Basis Cordes [®] RK	13A103	5
Ichthyol/Alliance healt...	Basis Cordes [®] RK	13A103	1
Ichthyol-Gesellschaft/V...	Basis Cordes [®] RK	13A103	1
Fagron	Basis Cordes [®] RK	13A103	2
Ichthyol-Gesellschaft C...	Basis Cordes [®] RK	13A103	1
Ichthyol-Gesellschaft /...	Basis Cordes [®] RK	13A103	1
Bombastus	Basis Cordes [®] RK	13A103	1
Ichthyol-Gesellschaft	Basis Cordes [®] RK	13A103	1
Sanacorp	Basis Cordes [®] RK	13A103	1
Ichthyol	Basis Cordes [®] RK	13A104	6
Ichthyol/Sanacorp	Basis Cordes [®] RK	13A104	3
Ichthyol/Gehe	Basis Cordes [®] RK	13A104	4
Ichtyol/Gehe	Basis Cordes [®] RK	13A206	1
Caelo	Basis Cordes [®] RK	13A206	2
Anzag/Ichtyol-Gesellsch...	Basis Cordes [®] RK	13A104	1
Caelo	Basis Cordes [®] RK	13A104	1
Ichthyol-Gesellschaft/N...	Basis Cordes [®] RK	13A206	2
Ichtyol-Gesellschaft	Basis Cordes [®] RK	13A206	1
Ichthyol/Sanacorp	Basis Cordes [®] RK	13A206	1
Ichthyol/Anzag	Basis Cordes [®] RK	13A206	1
Cordes, Hermann & Co	Basis Cordes [®] RK	13A206	2
Ichthyol-Gesellschaft/...	Basis Cordes [®] RK	13A206	1
Ichthyol-Gesellschaft /...	Basis Cordes [®] RK	13A206	1
Cordes, Hermani / Gehe	Basis Cordes [®] RK	13A206	1
Ichthyol	Basis Cordes [®] RK	13A206	4
AHD17.11.14 EK 16,06	Basis Cordes [®] RK	13A206	1
Ichthyol über Phoenix	Basis Cordes [®] RK	13A206	1
Ichtyol Gehe	Basis Cordes [®] RK	13A206	1
Ichtyol /Phönix	Basis Cordes [®] RK	13A206	1
VDL;8.12.14;15,91EUR	Basis Cordes [®] RK	13A206	2
Ichthyol-Gesellschaft	Basis Cordes [®] RK	13A206	1
Ichtyol - Gesellschaft	Basis Cordes [®] RK	13A206	1
Ichthyol/ Phoenix	Basis Cordes [®] RK	13A206	1
200	Basis Cordes [®] RK	13A206	1
Noweda/Wepa	Basis Cordes [®] RK	13A206	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
IchthyolAlliance Health...	Basis Cordes® RK	13A206	1
Phönix Gotha	Basis Cordes® RK	13A2C6	1
Ichthyol	Basis Cordes® RK	14A038	2
Noweda/Ichthyol	Basis Cordes® RK	14A038	1
Ichthyol/Gehe	Basis Cordes® RK	14A038	2
Ichthyol-Gesellschaft C...	Basis Cordes® RK	14A038	1
Ichthyol/Phönix	Basis Cordes® RK	14A038	4
Ichthyol/Sanacorp	Basis Cordes® RK	14A038	3
Ichthyol-Gesellschaft /...	Basis Cordes® RK	14A038	1
Ichthyol-Gesellschaft	Basis Cordes® RK	14A087	2
Anzag/Ichthyol-Gesellsch...	Basis Cordes® RK	14A051	1
Ichthyol	Basis Cordes® RK	14a087	1
Ichthyol /Gehe	Basis Cordes® RK	14A087	1
ICHTHYOL-GESELLSCHAFT C...	Basis Cordes® RK	14A087	1
Ichthyol Gesellschaft C...	Basis Cordes® RK	14A087	2
Ichthyol/ Noweda	Basis Cordes® RK	14A087	1
Ichthyol-Gesellschaft	Basis Cordes® RK	14A087	1
Ichthyol/Phönix	Basis Cordes® RK	14A087	1
Noweda	Basis Cordes® RK	14A087	1
AHD 04.02.15	Basis Cordes® RK	14A087	1
Ichthyol/Otto	Basis Cordes® RK	14A087	1
Ichthyol/Alliance Healt...	Basis Cordes® RK	14A087	1
Gehe/Ichthyol	Basis Cordes® RK	14A087	1
ICHTHYOL/GEHE	Basis Cordes® RK	14A087	1
Ichthyol	Basis Cordes® RK	14A087	1
Ichthyol/Phoenix	Basis Cordes® RK	14A148	4
Ichthyol/Sanacorp	Basis Cordes® RK	14A148	4
Ichthyol Gesellschaft	Basis Cordes® RK	14A148	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	14A148	2
Ichthyol/Phönix	Basis Cordes® RK	14A148	1
Ichthyol-Gesellschaft	Basis Cordes® RK	14A148	1
Ichthyol/Phoenix	Basis Cordes® RK	14A148	1
Ichthyol/Phoenix	Basis Cordes® RK	14A205	1
Ichthyol-Gesellschaft	Basis Cordes® RK	14A205	2
Ichthyol Gesellschaft	Basis Cordes® RK	14A205	1
Ichthyol-Gesellschaft	Basis Cordes® RK	14A205	5
Ichthyol/VDL	Basis Cordes® RK	14A205	1
Ichthyol/phönix	Basis Cordes® RK	14A205	3
Ichthyol/AHD	Basis Cordes® RK	14A205	1
Ichthyol / Phoenix	Basis Cordes® RK	14A205	1
Ichthyol	Basis Cordes® RK	14A205	2
Ichthyol/Fiebig	Basis Cordes® RK	14A205	1
Phönix Gotha	Basis Cordes® RK	14A205	1
ichthyol Gesellschaft ...	Basis Cordes® RK	14A205	1
Caelo	Basis Cordes® RK	14A205	3
Ichthyol/Noweda	Basis Cordes® RK	14A205	1
ICHTHYOL-GESELLSCHAFT	Basis Cordes® RK	14A205	1
Ichthyol-Gesellschaft /...	Basis Cordes® RK	14A205	1
Ichthyol/Gehe	Basis Cordes® RK	14A205	1
Caelo	Basis Cordes® RK	14D87	1
Ichthyol/Phoenix	Basis Cordes® RK	15A051	7
Ichthyol/Sanacorp	Basis Cordes® RK	15A051	5
Ichthyol/Anzag	Basis Cordes® RK	15A051	1
Ichthyol-Gesellschaft	Basis Cordes® RK	15A051	1
Phönix Gotha	Basis Cordes® RK	15A051	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	15A051	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	15A051	2
Caelo	Basis Cordes® RK	15A057	2
Ichthyol	Basis Cordes® RK	15A057	3
Ichthyol-Gesellschaft	Basis Cordes® RK	15A057	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol Gesellschaft	Basis Cordes [®] RK	15A057	3
Ichthyol-Gesellschaft	Basis Cordes [®] RK	15A057	7
Ichthyol/GEHE	Basis Cordes [®] RK	15A057	1
Ichthyol/Phoenix	Basis Cordes [®] RK	15A057	1
Ichthyol/ Sanacorp	Basis Cordes [®] RK	15A057	1
Ichthyol	Basis Cordes [®] RK	15A057	1
ICHTHYOL / 30	Basis Cordes [®] RK	15A057	1
ICHTHYOL-GESELLSCHAFT C...	Basis Cordes [®] RK	15A057	1
ICHTHYOL GESELLSCHAFT Co...	Basis Cordes [®] RK	15A057	1
Ichthyol Ges.	Basis Cordes [®] RK	15A057	1
Ichthyol-Gesellschaft/A...	Basis Cordes [®] RK	15A057	1
Ichthyol- Gesellschaft ...	Basis Cordes [®] RK	15A057	1
Noweda	Basis Cordes [®] RK	15A057	1
Ichthyol / Noweda	Basis Cordes [®] RK	15A057	1
Ichthyol/ Phönix	Basis Cordes [®] RK	15A115	1
Ichthyol/Gehe	Basis Cordes [®] RK	15A115	3
Ichthyol-Gesellschaft	Basis Cordes [®] RK	15A115	2
Ichthyol Gesellschaft C...	Basis Cordes [®] RK	15A115	1
Ichthyol/Sanacorp	Basis Cordes [®] RK	15A115	2
Ichthyol Gesellschaft	Basis Cordes [®] RK	15A115	1
Caelo	Basis Cordes [®] RK	15A115	2
Ichthyol	Basis Cordes [®] RK	15A115	2
ICHTHYOL-GESELLSCHAFT C...	Basis Cordes [®] RK	15A115	1
Ichthyol/N	Basis Cordes [®] RK	15A115	1
Noweda 29.032016	Basis Cordes [®] RK	15A115	2
Ichthyol Gesellschaft	Basis Cordes [®] RK	15A115	2
Ichthyol/Noweda	Basis Cordes [®] RK	15A115	1
Ichthyol / Sanacorp	Basis Cordes [®] RK	15A115	1
Ichthyol Werke	Basis Cordes [®] RK	15A115	1
Ichthyol/ Anzag	Basis Cordes [®] RK	15A115	1
phönix	Basis Cordes [®] RK	15A115	1
Ichthyol-Gesellschaft C...	Basis Cordes [®] RK	15A172	4
Ichthyol/Phönix	Basis Cordes [®] RK	15A172	3
Ichthyol-Gesellschaft C...	Basis Cordes [®] RK	15A172	4
Ichthyol	Basis Cordes [®] RK	15A172	1
Caelo	Basis Cordes [®] RK	15A172	1
Ichthyol/ Krieger	Basis Cordes [®] RK	15A172	1
Ichthyol-Gesellschaft	Basis Cordes [®] RK	15A172	1
Ichthyol/Sanacorp	Basis Cordes [®] RK	15A172	1
Ichthyol Gesellsch./Now...	Basis Cordes [®] RK	15A194	1
Ichthyol/Phönix	Basis Cordes [®] RK	15A194	4
Ichthyol-Gesellschaft C...	Basis Cordes [®] RK	15A194	1
Ichthyol/Noweda	Basis Cordes [®] RK	15A194	1
Ichthyol/Noweda	Basis Cordes [®] RK	15A194	1
Ichthyol / G	Basis Cordes [®] RK	15A194	1
Ichthyol/Sanacorp	Basis Cordes [®] RK	15A194	2
Caelo	Basis Cordes [®] RK	15A194	3
Noweda	Basis Cordes [®] RK	15A194	1
Ichthyol/Noweda	Basis Cordes [®] RK	15A194	1
Phönix	Basis Cordes [®] RK	15A194	2
Ichthyol/Alliance Healt...	Basis Cordes [®] RK	15A194	1
Ichthyol-Gesellschaft	Basis Cordes [®] RK	15A194	4
Euro OTC	Basis Cordes [®] RK	15A194	2
Ichthyol Werke	Basis Cordes [®] RK	15A194	1
ICHTHYOL / 20	Basis Cordes [®] RK	15A194	1
Ichthyol-Gesellschaft/AHD	Basis Cordes [®] RK	15A194	1
Ichthyol Gesellschaft/AHD	Basis Cordes [®] RK	15A194	1
Ichthyol	Basis Cordes [®] RK	15A057	1
ichthyol/noweda	Basis Cordes [®] RK	15A057	1
Ichthyol/Phoenix	Basis Cordes [®] RK	16A034	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol / Phönix	Basis Cordes® RK	16A034	1
Ichthyol-Gesellschaft /...	Basis Cordes® RK	16A034	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	16A034	1
Ichthyol	Basis Cordes® RK	16A034	4
Ichthyol/Phönix	Basis Cordes® RK	16A034	1
Cordes/Gehe	Basis Cordes® RK	16A034	1
Sanacorp	Basis Cordes® RK	16A034	1
Ichthyol-Gesellschaft	Basis Cordes® RK	16A034	1
Ichthyol/Noweda	Basis Cordes® RK	16A034	1
Ichthyol/Phönix	Basis Cordes® RK	16A034	1
Ichthyol/Sanacorp	Basis Cordes® RK	16A034	1
apomix/Noweda	Basis Cordes® RK	16A034	1
ICHTHYOL-GESELLSCHAFT	Basis Cordes® RK	16A035	1
Ichthyol-Gesellschaft	Basis Cordes® RK	16A035	2
Bombastus	Basis Cordes® RK	16A035	4
Ichthyol/Noweda	Basis Cordes® RK	16A035	1
Gehe/Ichthyol 23.11.16	Basis Cordes® RK	16A035	1
AHD/ichthyol	Basis Cordes® RK	16A035	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	16A035	1
Ichthyol	Basis Cordes® RK	16A035	1
Phoe	Basis Cordes® RK	16A035	1
Ichthyol	Basis Cordes® RK	16A035	2
Caelo	Basis Cordes® RK	16A035	1
Ichthyol/ Phoenix	Basis Cordes® RK	16A035	1
ICHTHYOL-GESELLSCHAFT/P...	Basis Cordes® RK	16A035	1
Ichthyol Gesellschaft	Basis Cordes® RK	16A035	1
Ichthyol-Gesellschaft/P...	Basis Cordes® RK	16A035	1
AHD 13.12.16 EK 16,54	Basis Cordes® RK	16A035	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	16A035	1
Gehe/Ichthyol	Basis Cordes® RK	16A035	1
Ichthyol/AHD	Basis Cordes® RK	16A035	1
Ichthyol Gesellschaft	Basis Cordes® RK	16A035	1
Ichthyol-Gesellschaft C...	Basis Cordes® RK	16A125	1
Ichthyol/Phönix	Basis Cordes® RK	16A125	5
ICHTHYOL GESELLSCHAFT Co...	Basis Cordes® RK	16A197	1
ichthyol/noweda	Basis Cordes® RK	17A046	4
Ichthyol-Gesellschaft	Basis Cordes® RK	17A046	2
Ichthyol/Noweda	Basis Cordes® RK	17A046	1
Ichthyol-Gesellschaft /...	Basis Cordes® RK	17A046	2
Ichthyol	Basis Cordes® RK	17A046	1
Ichthyol/Gehe	Basis Cordes® RK	17A046	1
Ichthyol/Noweda	Basis Cordes® RK	17A046	1
Ichthyol Gesellschaft/N...	Basis Cordes® RK	17A046	1
Ichthyol/ Noweda	Basis Cordes® RK	17A046	1
GEHE	Basis Cordes® RK	17A083	1
Ichthyol	Basis Cordes® RK	17A083	1
Ichthyol-Gesellschaft	Basis Cordes® RK	17A083	1
ichthyol/noweda	Basis Cordes® RK	17A083	1
Ichthyol	Basis Cordes® RK	17A149	1
ichthyol/noweda	Basis Cordes® RK	17A149	1
Ichthyol	Basis Cordes® RK	20061501	1
Ichthyol Gesellschaft/K...	Basis Cordes® RK	17A149	1
Ichthyol	Basis Cordes® RK	27102014B	1
Ichthyol	Basis Cordes® RK	28011507	1
Caelo	Basis Cordes® RK	P160/14	1
Ichthyol	Basis Cordes® RK	10A064	1
Ichthyol	Basis Cordes® RK	10A138	1
Ichthyol	Basis Cordes® RK	11A067	4
ICHTHYOL-GESELLSCHAFT	Basis Cordes® RK	11A079	1
ICHTHYOL-GES/GEHE	Basis Cordes® RK	11A079	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ebert	Basis Cordes [®] RK	11A079	1
Phönix Gotha	Basis Cordes [®] RK	11A079	1
Ichthyol	Basis Cordes [®] RK	11A079	1
Ichthyol Gesellschaft C...	Basis Cordes [®] RK	11A079	1
Ichthyol/ Noweda	Basis Cordes [®] RK	11A079	1
Ichthyol	Basis Cordes [®] RK	11A156	1
Ichthyol / Anzag	Basis Cordes [®] RK	11A156	1
Ichthyol Gesellschaft C...	Basis Cordes [®] RK	11A156	1
Noweda	Basis Cordes [®] RK	11A156	1
Ichtyol	Basis Cordes [®] RK	11A156	1
ichthyol	Basis Cordes [®] RK	1213Q-03300	1
Ichthyol - Gesellschaft	Basis Cordes [®] RK	11A079	1
Caelo	Basis Cordes [®] RK	12345	1
Ichthyol	Basis Cordes [®] RK	12A030	10
Ichthyol-Gesellschaft	Basis Cordes [®] RK	12A030	1
Gehe	Basis Cordes [®] RK	12A030	2
Ichthyol Gesellschaft C...	Basis Cordes [®] RK	12A045	1
Gehe	Basis Cordes [®] RK	12A045	4
Ichthyol-Gesellschaft C...	Basis Cordes [®] RK	12A045	1
Noweda	Basis Cordes [®] RK	12A045	1
Phönix Gotha	Basis Cordes [®] RK	12A045	1
Sanacorp/ Ichthyol Gese...	Basis Cordes [®] RK	12A045	1
Ichthyol	Basis Cordes [®] RK	12A045	2
Ichthyol-Gesellschaft / E	Basis Cordes [®] RK	12A045	1
Ichthyol Gesellschaft	Basis Cordes [®] RK	12A045	1
Ichthyol/Anzag	Basis Cordes [®] RK	12A045	1
ICHTHYOL-GES/GEHE	Basis Cordes [®] RK	12A045	1
Caelo	Basis Cordes [®] RK	12A045	1
Ichtyol 004	Basis Cordes [®] RK	12A046	1
Phönix/Ichthyol	Basis Cordes [®] RK	12A046	1
Ichthyol-Gesellschaft	Basis Cordes [®] RK	12A046	1
sanacorp	Basis Cordes [®] RK	12A046	1
Ichthyol/ Noweda	Basis Cordes [®] RK	12A046	1
Ichthyol Gesellschaft C...	Basis Cordes [®] RK	12A046	1
Ichthol	Basis Cordes [®] RK	12A046	1
Ichthyol	Basis Cordes [®] RK	12A046	1
Anzag/Ichtyol	Basis Cordes [®] RK	12A046	1
Noweda	Basis Cordes [®] RK	12A046	1
Ichthyol-Gesellschaft C...	Basis Cordes [®] RK	12A046	1
Ichthyol - Gesellschaft	Basis Cordes [®] RK	12A046	1
Anzag/ Ichtyol	Basis Cordes [®] RK	12A046	2
ICHTHYOL-GES.	Basis Cordes [®] RK	12A114	1
ICHTHYOL-GESELLSCHAFT	Basis Cordes [®] RK	12A046	1
Ichthyol	Basis Cordes [®] RK	12a046	1
Ichthyol-Gesellschaft	Basis Cordes [®] RK	12A114	1
Cordes/Phoenix	Basis Cordes [®] RK	12A114	1
Ichthyol	Basis Cordes [®] RK	12A114	5
Ichtyol	Basis Cordes [®] RK	12A114	1
Ichthyol	Basis Cordes [®] RK	12a114	1
Ichtyol-Gesellschaft/Ph...	Basis Cordes [®] RK	12A114	1
Phönix Gotha	Basis Cordes [®] RK	12A114	1
Ichthyol Werke	Basis Cordes [®] RK	12A114	1
Ichthyol/Ebert	Basis Cordes [®] RK	12A114	1
Ichthyol	Basis Cordes [®] RK	12a167	1
Ichthyol/Phönix	Basis Cordes [®] RK	12A114	1
Ichthyol	Basis Cordes [®] RK	12A167	6
Phönix	Basis Cordes [®] RK	12A167	2
Ichthyol/N/17.04.2013/P...	Basis Cordes [®] RK	12A167	1
Sanacorp / Ichthyol	Basis Cordes [®] RK	12A167	2
Ichthyol	Basis Cordes [®] RK	12ad30	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Cordes/Anzag	Basis Cordes® RK	13A038	2
Cordes/AHD	Basis Cordes® RK	13A038	1
Ichthyol	Basis Cordes® RK	13A038	9
Caelo	Basis Cordes® RK	13A038	2
Noweda	Basis Cordes® RK	13A038	1
Ichthyol/Sanacorp	Basis Cordes® RK	13A038	3
Ichthyol , Noweda	Basis Cordes® RK	13A038	1
Ichthyol -Gesellschaft	Basis Cordes® RK	13A038	1
Ichthyol-Gesellschaft	Basis Cordes® RK	13A038	1
Ichthyol GmbH, Phönix	Basis Cordes® RK	13A038	1
Ichthyol-Gesellschaft	Basis Cordes® RK	13A038	1
Alliance healthcare	Basis Cordes® RK	13A038	1
Ichtyol	Basis Cordes® RK	13A038	1
Bombastus	Basis Cordes® RK	13A038	1
Ichthyol-Gesellschaft / G	Basis Cordes® RK	13A038	1
Ichthyol GmbH	Basis Cordes® RK	13A038	1
Ichthyol/gehe	Basis Cordes® RK	13A039	1
Caelo	Basis Cordes® RK	13A039	2
Gehe	Basis Cordes® RK	13A039	1
Ichthyol, Noweda	Basis Cordes® RK	13A039	1
Ichthyol Gesellschaft	Basis Cordes® RK	13A039	2
Ichtyol /Jenne 16.01.2014	Basis Cordes® RK	13A039	1
ichthyol/AHD	Basis Cordes® RK	13A039	1
Ichthyol- Gesellschaft	Basis Cordes® RK	13A039	2
Gehe - Ichthyol	Basis Cordes® RK	13A039	1
Phoenix	Basis Cordes® RK	13A039	1
Ichthyol/gehe	Basis Cordes® RK	13A103	1
ichthyol	Basis Cordes® RK	13A103	1
Ichthyol/Noweda	Basis Cordes® RK	13A104	1
Ichthyol/N/PZN:01114881	Basis Cordes® RK	13A104	1
Ichthyol Gesellschaft	Basis Cordes® RK	13A104	1
Ichthyol Cordes	Basis Cordes® RK	13A104	1
Ichthyol/Gehe	Basis Cordes® RK	13A206	1
Ichthyol	Basis Cordes® RK	13A039	2
Ichthyol Gesellschaft C...	Basis Cordes® RK	15A057	1
ichthyol gesellschaft	Basis Cordes® RK	14A205	1
Ichthyol	Basis Cordes® RK	30041301	1
Ichthyol	Basis Cordes® RK	19111208	1
Caelo	Basis Cordes® RK	4312A-03300	1
Sonnenberg-Apotheke	Basis Cordes® RK	4643491	1

- 90 853 spectra from 501 *Apo-Ident* customers from a total of 26 955 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Basis Cordes® RK* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Basis Cordes® RK* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	140	0	35 404
Type B	0	130	0	93 579
Type C	0	477	22	90 579

The substance/substance group *Basis Cordes*[®] *RK* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.7143 %)
Type B	100.0000 % (> 99.9765 %)	100.0000 % (> 95.3846 %)
Type C	100.0000 % (> 99.6591 %)	95.5912 % (> 94.9900 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30892	30892	0.00	42.03
31078	30892	6.17	47.31
31660	31660	0.00	40.81

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Benzalkonium chloride solution 50%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30667-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Benzalkonium chloride solution 50%; Benzalkonii chloridi solutio 50%; Solutio benzalkonii chlorati 50%

Special notes

When selecting the *Benzalkonium chloride solution 50%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Benzalkonium chloride solution 50%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Benzalkonium chl...	11271919	30667	40	1402369
Euro OTC	Benzalkonium chl...	1105041V-03	31613	40	20130809
Caelo	Benzalkonium chl...	15051325	33605	40	20150804

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Benzalkonium chloride solution 50%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Benzalkonium chloride solution 50%*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Benzalkonium chloride soluti...	112719106	30543	40
Euro OTC	Benzalkonium chloride soluti...	1105041-01	30646	40
Euro OTC	Benzalkonium chloride soluti...	1105041V-03	31613 [†]	20
Caelo	Benzalkonium chloride soluti...	15051305	32960	40
Caelo	Benzalkonium chloride soluti...	15051327	33834	40

- 93 694 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 128 spectra from 90 *Apo-Ident* customers from 38 batches from the substance/substance group *Benzalkonium chloride solution 50%*.
- These include spectra of independent samples from 35 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Benzalkonium chloride soluti...	8041406	1
Euro OTC	Benzalkonium chloride soluti...	1105041V	1
Caelo	Benzalkonium chloride soluti...	13152105	6
Caelo	Benzalkonium chloride soluti...	13152111	4
Caelo	Benzalkonium chloride soluti...	13152114	1
Caelo	Benzalkonium chloride soluti...	13152115	1
Caelo	Benzalkonium chloride soluti...	13152123	3
Caelo	Benzalkonium chloride soluti...	13152121	1
Gehe	Benzalkonium chloride soluti...	13152123	1
Caelo	Benzalkonium chloride soluti...	13152125	18
Phönix	Benzalkonium chloride soluti...	13152125	1
Sanacorp	Benzalkonium chloride soluti...	13152125	1
Phönix Gotha	Benzalkonium chloride soluti...	13152125	1
Caelo	Benzalkonium chloride soluti...	15051305	8
Phoenix	Benzalkonium chloride soluti...	15051309	1
Caelo	Benzalkonium chloride soluti...	15051309	7
Caelo	Benzalkonium chloride soluti...	15051312	3
Caelo	Benzalkonium chloride soluti...	15051317	7
Caelo	Benzalkonium chloride soluti...	15051320	1
Caelo	Benzalkonium chloride soluti...	15051323	4
Caelo	Benzalkonium chloride soluti...	15051325	7
Ebert-Jacobi	Benzalkonium chloride soluti...	15051325	1
Caelo	Benzalkonium chloride soluti...	15119501	3
Fagron	Benzalkonium chloride soluti...	15119501	1
Caelo	Benzalkonium chloride soluti...	15119505	9
Caelo	Benzalkonium chloride soluti...	17044901	1
Sanacorp	Benzalkonium chloride soluti...	4/AKD2	1
Caesar/Noweda	Benzalkonium chloride soluti...	5012m-01168	1
Fagron	Benzalkonium chloride soluti...	7F21-N09	1
Caelo	Benzalkonium chloride soluti...	611benz	1
Fagron	Benzalkonium chloride soluti...	1105041-01	2
Euro OTC	Benzalkonium chloride soluti...	110041V-03	1
Euro OTC	Benzalkonium chloride soluti...	1105041-01	2
Gehe	Benzalkonium chloride soluti...	1105041-01	1
Euro OTC	Benzalkonium chloride soluti...	1105041V-03	2
Caelo	Benzalkonium chloride soluti...	11271906	2
Caelo	Benzalkonium chloride soluti...	11271912	1
Caelo	Benzalkonium chloride soluti...	11271919	4
Caelo	Benzalkonium chloride soluti...	12254603	1
Sanacorp	Benzalkonium chloride soluti...	12254608	1
Caelo	Benzalkonium chloride soluti...	12254610	2
Caelo	Benzalkonium chloride soluti...	12254608	2
Caelo	Benzalkonium chloride soluti...	12254615	2
Caelo	Benzalkonium chloride soluti...	12254620	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Benzalkonium chloride soluti...	12254615	1
Caelo	Benzalkonium chloride soluti...	1315225	1
Löhr Center/	Benzalkonium chloride soluti...	2013U-01168	1
Fagron	Benzalkonium chloride soluti...	2913A-01168	1
Caelo	Benzalkonium chloride soluti...	31011402	1

- 91 224 spectra from 501 *Apo-Ident* customers from a total of 26 973 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Benzalkonium chloride solution 50%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Benzalkonium chloride solution 50%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	180	0	93 694
Type C	0	127	1	91 224

The substance/substance group *Benzalkonium chloride solution 50%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.6667 %)
Type C	100.0000 % (> 99.6586 %)	99.2188 % (> 96.8750 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30667	30667	0.00	57.99
31613	30667	7.24	57.43
33605	30667	12.35	57.50

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Benzyl alcohol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30574-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Benzyl alcohol; Alkohol benzylicus

Special notes

When selecting the *Benzyl alcohol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Benzyl alcohol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Benzyl alcohol	13292103	31220	40	1403478
Caelo	Benzyl alcohol	13292107	31369	60	20131028
Caelo	Benzyl alcohol	16277302	32969	40	20161024

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 140 spectra of 3 reference samples from the substance/substance group *Benzyl alcohol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 404 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 130 spectra from 4 reference samples of the substance/substance group *Benzyl alcohol*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Benzyl alcohol	12108503	30574	40
Caelo	Benzyl alcohol	13292107	31369 [†]	10
Caelo	Benzyl alcohol	16341802	32968	40
Caelo	Benzyl alcohol	17191002	33857	40

- 93 744 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 29 spectra from 11 *Apo-Ident* customers from 16 batches from the substance/substance group *Benzyl alcohol*.
- These include spectra of independent samples from 14 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Benzyl alcohol	13292103	3
Caelo	Benzyl alcohol	13292107	2
Caelo	Benzyl alcohol	14206201	2
Caelo	Benzyl alcohol	15014901	4
Caelo	Benzyl alcohol	14206209	1
Sanacorp	Benzyl alcohol	15014907	1
Caelo	Benzyl alcohol	15014907	1
Caelo	Benzyl alcohol	15318403	3
Gehe	Benzyl alcohol	15318406	1
Caelo	Benzyl alcohol	15437503	3
Gehe	Benzyl alcohol	15437503	1
Gehe	Benzyl alcohol	15437509	1
Caelo	Benzyl alcohol	16341802	1
Gehe	Benzyl alcohol	15437512	1
Caelo	Benzyl alcohol	12108503	1
Caelo	Benzyl alcohol	12256208	1
Caelo	Benzyl alcohol	16341805	1
Caelo	Benzyl alcohol	13056803	1

- 91 323 spectra from 501 *Apo-Ident* customers from a total of 26 993 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Benzyl alcohol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Benzyl alcohol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	140	0	35 404
Type B	0	130	0	93 744
Type C	0	29	0	91 323

The substance/substance group *Benzyl alcohol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.7143 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.3846 %)
Type C	100.0000 % (> 99.6589 %)	100.0000 % (> 79.3103 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31220	31220	0.00	115.74
31369	31220	2.72	115.87
32969	31220	13.05	113.39

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Benzyl benzoate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31208-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Benzyl benzoate; Benzylum benzoicum

Special notes

When selecting the *Benzyl benzoate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Benzyl benzoate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Benzyl benzoate	13099320	31208	40	1403156
Caelo	Benzyl benzoate	14029602	31435	40	20140207
Caelo	Benzyl benzoate	171562	33695	40	20170530*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Benzyl benzoate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Benzyl benzoate*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Benzyl benzoate	09362308	30320	40
Caelo	Benzyl benzoate	14029602	31435 [†]	20
Caelo	Benzyl benzoate	16140805	32959	40
Caelo	Benzyl benzoate	172351	33696	40
Caelo	Benzyl benzoate	181819	34219	40

- 93 694 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 14 spectra from 7 *Apo-Ident* customers from 11 batches from the substance/substance group *Benzyl benzoate*.
- These include spectra of independent samples from 11 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caesar & Loretz GmbH	Benzyl benzoate	13099313	1
Caesar & Loretz GmbH	Benzyl benzoate	14029613	1
Caelo	Benzyl benzoate	15133401	1
Caesar & Loretz GmbH	Benzyl benzoate	14029603	1
Caesar & Loretz GmbH	Benzyl benzoate	15133401	1
Caelo	Benzyl benzoate	15243005	1
Bombastus	Benzyl benzoate	16140807	1
Caelo	Benzyl benzoate	11340415	1
Sanacorp	Benzyl benzoate	15243014	1
Caesar & Loretz GmbH	Benzyl benzoate	12141003	2
Phönix	Benzyl benzoate	12141008	1
Gehe	Benzyl benzoate	13099317	1
Caelo	Benzyl benzoate	14029613	1

- 91 338 spectra from 501 *Apo-Ident* customers from a total of 26 998 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Benzyl benzoate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Benzyl benzoate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	177	3	93 694
Type C	0	13	1	91 338

The substance/substance group *Benzyl benzoate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	98.3333 % (> 96.6667 %)
Type C	100.0000 % (> 99.6592 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31208	31208	0.00	44.64
31435	31208	2.13	47.12
33695	33695	0.00	41.28

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Benzyl nicotinate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30502-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Benzyl nicotinate; Benzyliis nicotinas; Benzylium nicotinicum

Special notes

When selecting the *Benzyl nicotinate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Benzyl nicotinate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Benzyl nicotinate	13272502	31199	40	1403109
Caelo	Benzyl nicotinate	14112404	31383	40	20140407
Caelo	Benzyl nicotinate	16358301	32967	40	20170113
Caelo	Benzyl nicotinate	163583	33055	40	20170113*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Benzyl nicotinate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 200 spectra from 5 reference samples of the substance/substance group *Benzyl nicotinate*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Benzyl nicotinate	12D25-N07	30502	40
Caelo	Benzyl nicotinate	11358501	30527	40
Caelo	Benzyl nicotinate	11358505	30528	40
Caelo	Benzyl nicotinate	15444804	32592	40
Caelo	Benzyl nicotinate	16358302	32966	40

- 93 674 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 47 spectra from 31 *Apo-Ident* customers from 25 batches from the substance/substance group *Benzyl nicotinate*.
- These include spectra of independent samples from 23 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Beiersdorf AG	Benzyl nicotinate	15293104	1
Anzag	Benzyl nicotinate	12282802	1
Caelo	Benzyl nicotinate	10111403	1
Caelo	Benzyl nicotinate	13157201	1
Caelo	Benzyl nicotinate	13272502	2
Caelo	Benzyl nicotinate	14112408	4
Caelo	Benzyl nicotinate	14112405	1
Caelo	Benzyl nicotinate	12142008	2
Caelo	Benzyl nicotinate	14112404	1
Caelo	Benzyl nicotinate	15124402	4
Caelo	Benzyl nicotinate	1529304	1
Caelo	Benzyl nicotinate	15124408	1
Caelo	Benzyl nicotinate	15293104	1
Caelo	Benzyl nicotinate	15444803	1
Caelo	Benzyl nicotinate	15444804	3
Caelo	Benzyl nicotinate	16358302	2
Caelo	Benzyl nicotinate	16154409	1
Caelo	Benzyl nicotinate	16358308	3
Caelo	Benzyl nicotinate	2091301	1
Anzag 10.08.2012	Benzyl nicotinate	12142004	1
Caelo	Benzyl nicotinate	12142004	1
anzag	Benzyl nicotinate	12282802	1
Caelo	Benzyl nicotinate	12282802	2
VDL	Benzyl nicotinate	12D25-N07	1
Caelo	Benzyl nicotinate	12282804	1
Fagron	Benzyl nicotinate	12I20-N01	1
Anzag/ichthyol	Benzyl nicotinate	13035408	1
Caelo	Benzyl nicotinate	13035402	2
Caelo	Benzyl nicotinate	13157205	3
Phönix 28.01.2014	Benzyl nicotinate	13272502	1

- 91 305 spectra from 501 *Apo-Ident* customers from a total of 26 984 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Benzyl nicotinate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all

relevant spectra of the various substances were compared with *Benzyl nicotinate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	200	0	93 674
Type C	0	46	1	91 305

The substance/substance group *Benzyl nicotinate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.0000 %)
Type C	100.0000 % (> 99.6587 %)	97.8723 % (> 91.4894 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31199	31199	0.00	32.06
31383	31199	3.25	29.76
32967	33055	1.63	33.54
33055	33055	0.00	33.82

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Bergamot oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31678-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Bergamot oil; Oleum bergamottae; Oleum citrus bergamia

Special notes

When selecting the *Bergamot oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Bergamot oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Bergamot oil	15042001	31678	40	AR-15-FG-008850-02
Caelo	Bergamot oil	15042001	31696	40	20150217

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Bergamot oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 1 different batches.
- 35 464 spectra from a total of 839 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 212 spectra from 7 reference samples of the substance/substance group *Bergamot oil*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Bergamot oil	15042001	31678 [†]	20
Caelo	Bergamot oil	15042001	31696 [†]	20
Caelo	Bergamot oil	16083611	33497	40
Caelo	Bergamot oil	16083611	33655	40
Caelo	Bergamot oil	171992	33894	40
Caelo	Bergamot oil	180207	34230	40
Caelo	Bergamot oil	16083601	not required	12

- 93 662 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 7 spectra from 7 *Apo-Ident* customers from 6 batches from the substance/substance group *Bergamot oil*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Bergamot oil	15042001	2
Caelo	Bergamot oil	15042008	1
Bombastus	Bergamot oil	300470	1
Caelo	Bergamot oil	16083601	1
Taoasis	Bergamot oil	19821-122171	1
Taoasis	Bergamot oil	53731-127819	1

- 91 345 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Bergamot oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Bergamot oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	206	6	93 662
Type C	0	7	0	91 345

The substance/substance group *Bergamot oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	97.1698 % (> 95.7547 %)
Type C	100.0000 % (> 99.6598 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several

new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31678	31678	0.00	60.23
31696	31678	1.34	60.98

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Cajeput
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31502-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Cajeput; Oleum Cajeputi artificiale

Special notes

When selecting the *Cajeput* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cajeput*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cajeput	13141409	31502	40	1412341
Caelo	Cajeput	14225206	31605	40	20140730

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Cajeput*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 160 spectra from 5 reference samples of the substance/substance group *Cajeput*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cajeput	13141409	31502 [†]	20
Caelo	Cajeput	14225206	31605 [†]	20
Caelo	Cajeput	162770	32688	40
Caelo	Cajeput	16277003	33620	40
Caelo	Cajeput	172410	33910	40

- 93 714 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 8 spectra from 3 *Apo-Ident* customers from 7 batches from the substance/substance group *Cajeput*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caesar & Loretz/Noweda	Cajeput	14225210	1
Caelo	Cajeput	15063501	2
Fagron	Cajeput	15063511	1
Taoasis	Cajeput	1518-122595	1
Caelo	Cajeput	16277004	1
Taoasis	Cajeput	6781012-119058	1
Taoasis	Cajeput	6781010-120111	1

- 91 344 spectra from 501 *Apo-Ident* customers from a total of 27 002 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cajeput* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cajeput* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	160	0	93 714
Type C	0	7	1	91 344

The substance/substance group *Cajeput* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.2500 %)
Type C	100.0000 % (> 99.6597 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31502	31502	0.00	21.39
31605	31502	11.23	19.63

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Camphor oil 20%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30551-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Camphor oil 20%; Oleum camphoratum 20%

Special notes

When selecting the *Camphor oil 20%* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Camphor oil 20%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Camphor oil 20%	12261202	31154	40	1404534
Caelo	Camphor oil 20%	13369902	31347	40	20131112

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Camphor oil 20%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Camphor oil 20%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Camphor oil 20%	11121003	30551	40
Caelo	Camphor oil 20%	13369902	31347 [†]	20
Caelo	Camphor oil 20%	17083502	33494	40
Caelo	Camphor oil 20%	17083502	33495	40

- 93 734 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 54 spectra from 27 *Apo-Ident* customers from 17 batches from the substance/substance group *Camphor oil 20%*.
- These include spectra of independent samples from 16 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Camphor oil 20%	13369902	4
Noweda	Camphor oil 20%	14147303	1
Caelo	Camphor oil 20%	15246805	2
Caelo	Camphor oil 20%	15246809	6
Caelo	Camphor oil 20%	16316101	2
Caelo	Camphor oil 20%	17301902	1
Phoenix	Camphor oil 20%	16316101	1
Caelo	Camphor oil 20%	17083502	1
Caelo	Camphor oil 20%	9250001	3
Caelo	Camphor oil 20%	12122502	9
ahd/lamotte	Camphor oil 20%	133699902	1
Caelo	Camphor oil 20%	14028202	2
Caelo	Camphor oil 20%	15056907	11
Caelo	Camphor oil 20%	15056906	1
gehe	Camphor oil 20%	15246805	1
Caelo	Camphor oil 20%	16085301	1
Caelo	Camphor oil 20%	12261203	2
Phönix	Camphor oil 20%	13039603	1
Caelo	Camphor oil 20%	13039603	1
Sanacorp	Camphor oil 20%	13039603	1
Caelo	Camphor oil 20%	1226103	1
Caelo	Camphor oil 20%	14147303	1

- 91 298 spectra from 501 *Apo-Ident* customers from a total of 26 992 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Camphor oil 20%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Camphor oil 20%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	140	0	93 734
Type C	0	53	1	91 298

The substance/substance group *Camphor oil 20%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6587 %)	98.1481 % (> 92.5926 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31154	31154	0.00	29.37
31347	31154	2.24	28.88

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Camphor ointment**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 32515-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Camphor ointment; Unguentum camphoratum

Special notes

When selecting the *Camphor ointment* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Camphor ointment*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Camphor ointment	161773	32515	40	20160711*
Caelo	Camphor ointment	17045702	33224	40	20170214

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Camphor ointment*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 40 spectra from 1 reference samples of the substance/substance group *Camphor ointment*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Camphor ointment	180955	34159	40

- 93 834 spectra from a total of 2166 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 0 spectra from 0 *Apo-Ident* customers from 0 batches from the substance/substance group *Camphor ointment*.
- These include spectra of independent samples from 0 batches, of which no spectra have been used for the generation of the database.
- 91 352 spectra from 501 *Apo-Ident* customers from a total of 27 009 batches from a further 334 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Camphor ointment* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Camphor ointment* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	40	0	93 834
Type C	0	0	0	91 352

The substance/substance group *Camphor ointment* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 85.0000 %)
Type C	100.0000 % (> 99.6586 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
32515	32515	0.00	21.73
33224	32515	7.87	26.87

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Camphor spirit
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30937-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Camphor spirit; Solutio Camphorae spirituosae; Spiritus camphoratus

Special notes

When selecting the *Camphor spirit* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Camphor spirit*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Camphor spirit	12196804	30937	50	20120721
Caelo	Camphor spirit	13354201	31389	40	1506005
Euro OTC	Camphor spirit	1506006	31700	40	20150625*
Euro OTC	Camphor spirit	1506007	31701	40	20150625*
Caelo	Camphor spirit	16039501	32307	40	20160229

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 210 spectra of 5 reference samples from the substance/substance group *Camphor spirit*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 334 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 4 reference samples of the substance/substance group *Camphor spirit*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Camphor spirit	10058403	30410	40
Euro OTC	Camphor spirit	1506006	31700 [†]	20
Euro OTC	Camphor spirit	1506007	31701 [†]	20
Caelo	Camphor spirit	170923	33327	40

- 93 754 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 79 spectra from 47 *Apo-Ident* customers from 39 batches from the substance/substance group *Camphor spirit*.
- These include spectra of independent samples from 36 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
GEHE	Camphor spirit	1301024-02	1
Caelo	Camphor spirit	10514C	1
Euro OTC	Camphor spirit	1307040	1
Caelo	Camphor spirit	13354204	2
Caelo	Camphor spirit	5031507	1
Caelo	Camphor spirit	13354206	1
Caelo	Camphor spirit	13354210	1
Euro OTC	Camphor spirit	1405026-01	1
Euro OTC	Camphor spirit	1405029-01	4
Eu Rho	Camphor spirit	1405029-01	1
Euro OTC	Camphor spirit	1405030	1
Euro OTC	Camphor spirit	140503010	1
Euro OTC	Camphor spirit	1405030-01	2
Euro OTC	Camphor spirit	1506006-01	5
Caelo	Camphor spirit	16039505	2
Caelo	Camphor spirit	16039501	1
Caelo	Camphor spirit	16141105	3
Caelo	Camphor spirit	17092302	1
gepepharm/Gehe	Camphor spirit	17092303	1
Caelo	Camphor spirit	17092303	1
Caelo	Camphor spirit		1
Caelo	Camphor spirit	17318502	1
Euro OTC	Camphor spirit	110022-02	1
Caelo	Camphor spirit	10058403	2
Caelo	Camphor spirit	11347401	3
Euro OTC	Camphor spirit	1108022-02	1
Anzag ek:13.14EUR 14.08.12	Camphor spirit	11347401	1
Sanacorp	Camphor spirit	1210026-01	1
Euro OCT	Camphor spirit	1210028-01	1
Euro OTC	Camphor spirit	1210028-01	6
anzag, 05.12.2012	Camphor spirit	11347401	1
ANZAG	Camphor spirit	121002801	1
Caelo	Camphor spirit	12196803	6
Euro OTC	Camphor spirit	121028-01	1
Spreewälder/AHD	Camphor spirit	1301024-02	1
Caelo	Camphor spirit	12196804	1
Euro OTC	Camphor spirit	1301024-02	2
Euro OTC	Camphor spirit	1307039-01	4
Phoenix 03.12.12	Camphor spirit	1307039-01	1
Euro OTC	Camphor spirit	1307040-01	1
Caelo	Camphor spirit	140503001	1
Caelo	Camphor spirit	22101201	1
Caelo	Camphor spirit	62143436	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Camphor spirit	30101309	1
Euro OTC	Camphor spirit	19071208	1
SprewälderArzneimittel ...	Camphor spirit	702110	1
PKH GmbH	Camphor spirit	1405030-01	1
Euro OTC	Camphor spirit	1108022-01	1
Caleo/ Anzag	Camphor spirit	11347401	1

- 91 273 spectra from 501 *Apo-Ident* customers from a total of 26 971 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Camphor spirit* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Camphor spirit* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	210	0	35 334
Type B	0	120	0	93 754
Type C	0	73	6	91 273

The substance/substance group *Camphor spirit* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.1429 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6587 %)	92.4051 % (> 88.6076 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30937	31700	5.47	95.22
31389	31389	0.00	101.46
31700	31700	0.00	96.38
31701	31701	0.00	98.86
32307	31389	5.84	101.98

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Caraway oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30876-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Caraway oil; Carvi aetheroleum; Oleum carvi

Special notes

When selecting the *Caraway oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Caraway oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Caraway oil	12J22-N03	30876	39	20121025
Caelo	Caraway oil	12043508	31024	40	AR-14-FG-012307-01
Caelo	Caraway oil	14166811	31666	40	20140604

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 119 spectra of 3 reference samples from the substance/substance group *Caraway oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 425 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Caraway oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Caraway oil	14166811	31666 [†]	20
Caelo	Caraway oil	16252906	33252	40
Caelo	Caraway oil	181649	34229	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 161 spectra from 84 *Apo-Ident* customers from 78 batches from the substance/substance group *Caraway oil*.
- These include spectra of independent samples from 75 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Caraway oil	14166803	3
Bombastus	Caraway oil	284758	4
Bombastus	Caraway oil	280523	1
Caelo	Caraway oil	2031504	1
Caelo	Caraway oil	5091402	1
Caesar/Noweda	Caraway oil	712Q-01968	2
Caelo	Caraway oil	10319108	2
Caelo	Caraway oil	10319115	2
Caelo	Caraway oil	6091401	1
Caelo	Caraway oil	10319116	1
Noweda	Caraway oil	10319109	1
Caelo	Caraway oil	10319117	1
Caelo	Caraway oil	10319118	5
Caelo	Caraway oil	10319120	2
Schubert	Caraway oil	1031E0-801	1
Fagron	Caraway oil	11/13-N03	2
Fagron	Caraway oil	11A19-N08	1
Fagron	Caraway oil	11L13-N03	2
Großhandel	Caraway oil	11L13N03	1
Caelo	Caraway oil	12013508	1
Caelo	Caraway oil	12043502	2
Caelo	Caraway oil	12043503	4
Caelo	Caraway oil	12043507	3
Caelo	Caraway oil	12043508	3
Caelo	Caraway oil	12043511	9
Phönix	Caraway oil	12043511	1
Caelo	Caraway oil	12043513	8
Caelo	Caraway oil	12043512	1
Caelo	Caraway oil	12043516	1
Fagron	Caraway oil	12J22-N03	1
Phönix	Caraway oil	12J22-N03	1
Caelo	Caraway oil	12043518	1
Fagron	Caraway oil	12J2-N03	1
Fagron	Caraway oil	13C08-N03	3
Fagron	Caraway oil	13c08-N03	1
Fagron	Caraway oil	13E06-N13	3
Fagron	Caraway oil	13c08-N08	1
Großhandel	Caraway oil	13E06-N13	1
Caelo	Caraway oil	14166802	4
Caelo	Caraway oil	14166806	2
Caelo	Caraway oil	14166810	2
Caelo	Caraway oil	14166811	2
Caelo	Caraway oil	14166812	1
Caelo	Caraway oil	1508506	1
Caelo	Caraway oil	15085703	1
Caelo	Caraway oil	15085704	3
Caelo	Caraway oil	15085706	2
Caelo	Caraway oil	15085707	4
Caelo	Caraway oil	15085709	2
Caelo	Caraway oil	15085710	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Caraway oil	15085713	1
Caelo	Caraway oil	16252902	1
Caelo	Caraway oil	16252903	3
Caelo	Caraway oil	16252906	2
Caelo	Caraway oil	17101301	1
Caelo	Caraway oil	22051306	1
Caelo	Caraway oil	19121207	1
Bombastus	Caraway oil	268955	1
Caelo	Caraway oil	27041503	1
Anzag 16.06.2010	Caraway oil	2410Q-01968	1
Bombastus	Caraway oil	270424	3
Klenk	Caraway oil	274361	1
Bombastus	Caraway oil	274361	2
Bombastus	Caraway oil	272320	1
Bombastus	Caraway oil	276856	3
Caelo	Caraway oil	28091205	1
Phönix	Caraway oil	276856	1
phönix 14.08.2013	Caraway oil	276856	1
Bombastus	Caraway oil	287725	1
Fagron	Caraway oil	287725	1
Bombastus	Caraway oil	290759	1
Bombastus	Caraway oil	291934	4
Bombastus	Caraway oil	289019	1
Bombastus	Caraway oil	295058	3
Bombastus	Caraway oil	297505	3
Bombastus	Caraway oil	295572	1
Caelo	Caraway oil	300797	1
Bombastus	Caraway oil	303305	2
Caelo	Caraway oil	297505	1
Bombastus	Caraway oil	300797	1
Bombastus	Caraway oil	308593	1
Fagron	Caraway oil	628194	1
Fagron	Caraway oil	699440	1
Fagron	Caraway oil	6999440	1
	Caraway oil	81816288	1
Bombastus	Caraway oil	299609	1
Bombastus	Caraway oil	287889	1

- 91 191 spectra from 501 *Apo-Ident* customers from a total of 26 932 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Caraway oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Caraway oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	119	0	35 425
Type B	0	100	0	93 774
Type C	0	157	4	91 191

The substance/substance group *Caraway oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 94.9580 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6586 %)	97.5155 % (> 95.6522 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30876	31024	7.55	32.83
31024	31024	0.00	36.14
31666	31024	11.18	33.59

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Castor oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30432-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Castor oil; Oleum ricini; Ricini oleum

Special notes

When selecting the *Castor oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of castor oil or rizol raw material: Castor oil and rizol (castor oil enriched with oxygen) can not be conclusively distinguished by spectral analysis. An additional analysis of identity may be required.

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Castor oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Castor oil	12114603	30741	40	20120504
Caelo	Castor oil	14020104	31343	40	AR-15-FG-005211-01
Euro OTC	Castor oil	1509003	31739	45	20150918*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Castor oil	14304807	31751	40	20141106
Fagron	Castor oil	15D02-B02-314240	31824	40	20150410
Euro OTC	Castor oil	1512024	31845	40	20160111*
Fagron	Castor oil	16J12-B01-328159	32943	40	20161121

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 285 spectra of 7 reference samples from the substance/substance group *Castor oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 7 different batches.
- 35 259 spectra from a total of 833 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 795 spectra from 22 reference samples of the substance/substance group *Castor oil*.
- These include spectra of independent samples from 18 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Bombastus	Castor oil	00266593	30432	40
Euro OTC	Castor oil	1109001-01	30608	40
Caelo	Castor oil	12149306	30776	40
Euro OTC	Castor oil	1509003	31739 [†]	15
Caelo	Castor oil	14304807	31751 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Castor oil	15D02-B02-314240	31824 [†]	20
Euro OTC	Castor oil	1512024	31845 [†]	20
Euro OTC	Castor oil	1604019	32192	40
Caelo	Castor oil	161028	32201	40
Caelo	Castor oil	160922	32207	40
Caelo	Castor oil	162000	32311	40
Caelo	Castor oil	161672	32312	40
Caelo	Castor oil	162362	32344	40
Caelo	Castor oil	162731	32681	40
Euro OTC	Castor oil	1611005	32858	40
Caelo	Castor oil	170272	33071	40
Caelo	Castor oil	163146	33072	40
Caelo	Castor oil	170271	33073	40
Euro OTC	Castor oil	1705035	33596	40
Caelo	Castor oil	171899	33678	40
Euro OTC	Castor oil	1801007	34004	40
Euro OTC	Castor oil	1807003-01	34211	40

- 93 079 spectra from a total of 2145 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1577 spectra from 313 *Apo-Ident* customers from 348 batches from the substance/substance group *Castor oil*.
- These include spectra of independent samples from 340 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Euro OTC	Castor oil	1407053-02	10
Caelo	Castor oil	11359502	1
GEHE	Castor oil	11359507	1
Caelo	Castor oil	12114603	2
Euro OTC	Castor oil	121008-01	1
Caelo	Castor oil	12114605	11
Caelo	Castor oil	12114606	2
Phönix	Castor oil	1212008-01	1
Euro OTC	Castor oil	1212008-01	4
Euro OTC	Castor oil	1212008-02	8
Caelo	Castor oil	12270106	8
Caelo	Castor oil	12270108	3
Caelo	Castor oil	12270111	5
Fagron	Castor oil	12H29-N04	2

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Castor oil	1304013-01	9
EuRho	Castor oil	1304013-01	1
Euro OTC	Castor oil	1308019-01	11
Euro OTC	Castor oil	1304013-02	11
Caelo	Castor oil	13096304	10
Caelo	Castor oil	13134604	22
Caelo	Castor oil	13134609	14
gehe 25.07.2013	Castor oil	275937	1
Fagron	Castor oil	13C08-N08	2
Caelo	Castor oil	278661	2
Euro OTC	Castor oil	3413U03553	1
Fargon /12.07.2013 / No...	Castor oil	13D29-N01	1
Henry Lamotte oils GmbH...	Castor oil	10000000714	1
Euro OTC	Castor oil	1407053-01	3
Euro OTC	Castor oil	1412006-02	11
Euro OTC	Castor oil	1412006-01	5
Caelo	Castor oil	14144305	2
Caelo	Castor oil	1504006-02	2
Caelo	Castor oil	14199603	4
Euro OTC	Castor oil	1509003-02	15
Euro OTC	Castor oil	1509003-01	21
Euro OTC	Castor oil	1512024	1
Euro OTC	Castor oil	1512024-02	9
Caelo	Castor oil	15144304	5
Caelo	Castor oil	15392001	11
Caelo	Castor oil	16016203	6
Caelo	Castor oil	16016206	4
Euro OTC	Castor oil	1604019-01	14
Caelo	Castor oil	16016207	4
Caelo	Castor oil	290212	4
Sanacorp;09,04.15; 3,33EUR	Castor oil	291722	1
Bombastus	Castor oil	292409	13
HenryLamotte/Anzag	Castor oil	2114I-02009	1
Bombastus	Castor oil	290212	6
VDL; 9.58EUR; 6.10.2015	Castor oil	294667	2
Bombastus	Castor oil	294668	14
Bombastus	Castor oil	295412	17
Bombastus	Castor oil	296751	12
Bombastus	Castor oil	298698	4
Bombastus	Castor oil	298699	23
Bombastus	Castor oil	300629	14
Bombastus	Castor oil	301590	3
Sanacorp, 9,48EUR, 17.1.17	Castor oil	301587	1
Caelo	Castor oil	305230	1
Jenne,Lamotte	Castor oil	123384	1
Caelo	Castor oil	13043601	13
Anzag	Castor oil	5334001	4
Bombastus	Castor oil	273390	5
Bombastus	Castor oil	270580	1
Phönix	Castor oil	123384	1
Lamotte / Gehe	Castor oil	123384	2
Noweda / Henry Lamotte	Castor oil	123384	2
Lamotte/Anzag	Castor oil	123384L902	1
Otto Fischar	Castor oil	1000001914	1
Lamotte	Castor oil	123812	2
Caelo	Castor oil	12121307	1
Euro OTC	Castor oil	1206039-01	4
Caelo	Castor oil	12149306	7
Caelo	Castor oil	12318204	1
Caelo	Castor oil	12318207	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Castor oil	13043601	1
Caelo	Castor oil	22071302	1
Bombastus	Castor oil	273882	9
Caelo	Castor oil	22051303	1
Bombastus	Castor oil	275937	10
Caelo	Castor oil	13065405	28
Bombastus	Castor oil	275938	12
Anzag 15.03.2013	Castor oil	275938	1
Bombastu-Werke AG/Gehe	Castor oil	278660	1
Anzag 06.05.2013	Castor oil	278661	2
Bombastus	Castor oil	278661	13
Bombastus	Castor oil	281943	11
Caelo	Castor oil	30031302	1
Bombastus	Castor oil	280791	17
Noweda / Henry Lamotte	Castor oil	5334001	2
Caelo	Castor oil	5334001	1
Henry Lamotte Oils/Gehe	Castor oil	5334001	2
Lamotte / Gehe	Castor oil	5334001	1
phö	Castor oil	5450801	1
Lamotte	Castor oil	5450801	4
Henry Lamotte/Ebert	Castor oil	5450801	3
Caelo	Castor oil	5450801	4
Henry Lamotte	Castor oil	5450801	1
Caelo	Castor oil	13196309	8
Euro OTC	Castor oil	1403024-01	9
Caelo	Castor oil	14058810	5
Caelo	Castor oil	14199501	13
Gutenberg-Apotheke	Castor oil	281944	2
Caelo	Castor oil	14199502	3
Bombastus	Castor oil	281944	6
VDL;17.10.14.3,96EUR	Castor oil	288639	1
Bombastus	Castor oil	291722	14
VDL	Castor oil	290213	1
Caelo	Castor oil	P161/14	1
Caelo	Castor oil	13134603-3417310103	1
Caelo	Castor oil	13134608-3438311151	1
Fagron	Castor oil	13134609	1
Caelo	Castor oil	13196313	13
Caelo	Castor oil	13G29-B01	1
Henry Lamotte/Sanacorp	Castor oil	123384	3
Lamotte/N/04.06.13/PZN:...	Castor oil	123384	1
Lamotte Phönix	Castor oil	5450801	2
Sanacorp	Castor oil	12270111	1
Caelo	Castor oil	1091505	1
Henry Lamotte / Noweda	Castor oil	123384	1
Caelo	Castor oil	123384	2
Euro OTC	Castor oil	1112013C	1
Henry Lamotte Oils GmbH...	Castor oil	123384	1
Noweda/Henry Lamotte	Castor oil	123812	4
Sanacorp	Castor oil	123812	3
Phoenix/Lamotte	Castor oil	123812	2
Henry Lamotte Oils/Noweda	Castor oil	123812	1
Henry Lamotte	Castor oil	123812	1
Henry Lamotte/Noweda	Castor oil	123812	1
Lamotte/Anzag	Castor oil	123812	1
Lamotte/Gehe	Castor oil	123812	2
H.Lamotte/Noweda	Castor oil	123812	1
Lamotte/Sanacorp	Castor oil	123812	1
Henry Lamotte Oils GmbH	Castor oil	123812	1
Stumpf/Henry Lamotte	Castor oil	124189	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Castor oil	124189	3
Lamotte/Sanacorp	Castor oil	124189	3
Henry Lamotte Oils GmbH...	Castor oil	124189	1
Lamotte	Castor oil	124189	2
lamotte/Gehe	Castor oil	124189	2
Henry Lamotte/Sanacorp	Castor oil	124189	1
Lamotte/Phoenix	Castor oil	124189	1
Henry Lamotte Oils/GEHE	Castor oil	124637	1
Henry Lamotte Oils GmbH.	Castor oil	124189	1
lamotte/Gehe	Castor oil	124637	1
Caelo	Castor oil	124637	1
Lamotte/Alliance	Castor oil	124637	2
Noweda/Henry Lamotte	Castor oil	124637	1
Henry Lamotte/Noweda	Castor oil	124637	1
Lamotte/Sanacorp	Castor oil	125162	1
H.Lamotte/Noweda	Castor oil	125162	2
Hedinger	Castor oil	125162	1
Stumpf/Henry Lamotte	Castor oil	125162	1
Phoenix/ Henry Lamotte	Castor oil	125162	1
Lamotte / Noweda	Castor oil	125162	1
Henry la Motte	Castor oil	125162	1
lamotte phönix	Castor oil	125162	1
Henry Lamotte/Phoenix	Castor oil	125473	2
Noweda/Henry Lamotte	Castor oil	125473	1
Lamotte/Sanacorp	Castor oil	125473	2
Euro OTC	Castor oil	2022015C	1
Bombastus	Castor oil	415Q-02009	1
Caelo	Castor oil	6051411	1
gehe	Castor oil	6861795	1
Bombastus	Castor oil	6861789	1
Henry Lamotte / Kehr	Castor oil	1000000714	1
Ebert-Jacobi	Castor oil	6861795	1
Henry Lamotte Oils GmbH	Castor oil	1000000714	1
Henry Lamotte	Castor oil	1000000714	1
Henry Lamotte/Phoenix	Castor oil	1000000714	1
Henry Lamotte/AHD	Castor oil	1000001128	2
Henry Lamotte	Castor oil	1000001128	3
Henry Lamotte oils GmbH...	Castor oil	1000001128	1
Henry Lamotte/Gehe	Castor oil	1000001128	1
LaMotte / Phoenix	Castor oil	1000001128	3
Lamotte/Sanacorp	Castor oil	1000001128	2
Lamotte/Anzag	Castor oil	1000001914	2
Noweda/Henry Lamotte	Castor oil	1000001914	1
Henry Lamotte	Castor oil	1000001914	4
H.Lamotte/Noweda	Castor oil	1000001914	1
Lamotte/Sanacorp	Castor oil	1000001914	3
Henry Lamotte/Phoenix	Castor oil	1000001914	3
Henry Lamotte/Noweda	Castor oil	1000001914	1
Henry Lamotte Oils GmbH	Castor oil	1000001914	1
Henry Lamotte/Noewda	Castor oil	100000238	1
Lamotte/Anzag	Castor oil	1000002384	1
Henry Lamotte Oils GmbH	Castor oil	1000002384	1
Lamotte/Sanacorp	Castor oil	1000002384	2
Caelo	Castor oil	1000002384	1
Henry Lamotte	Castor oil	1000002384	2
LaMotte / Phoenix	Castor oil	1000002384	1
Henry Lamotte/Noweda	Castor oil	1000002384	1
Alte Apotheke	Castor oil	1000002384	1
Henry Lamotte	Castor oil	1000002766	3
LaMotte / Phoenix	Castor oil	1000002766	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte/ Phönix	Castor oil	1000003564	1
Caelo	Castor oil	1000003564	5
Lamotte/Phönix	Castor oil	1000003564	4
Henry Lamotte	Castor oil	1000003564	4
Noweda	Castor oil	1000003564	1
lamotte/ Phoenix	Castor oil	1000003564	1
Lamotte/Sanacorp	Castor oil	1000004198	1
Henry Lamotte/ Phönix	Castor oil	1000004198	1
lamotte/ Phoenix	Castor oil	1000004198	1
Henry Lamotte	Castor oil	1000004198	1
Lamotte/Phönix	Castor oil	1000004198	1
lamotte	Castor oil	1000004198	1
Noweda / H. Lamotte	Castor oil	1000004198	1
Henry Lamotte/Noewda	Castor oil	1000004198	1
lamotte/ Phoenix	Castor oil	1000004554	1
H.Lamotte	Castor oil	1000004554	2
Fiebig	Castor oil	1000004554	1
Henry Lamotte/Noweda	Castor oil	1000004554	3
Lamotte/Phönix	Castor oil	1000005199	1
Lamotte/Noweda	Castor oil	1000005199	1
Lamotte/Phönix	Castor oil	1000005818	2
Henry la Motte	Castor oil	1000005818	1
Kehr/Henry Lamotte	Castor oil	1000005818	1
Caelo	Castor oil	1000005818	1
Caelo	Castor oil	100004198	1
Caelo	Castor oil	10021402	1
Euro OTC	Castor oil	10022015D	1
Euro OTC	Castor oil	1014A03553	1
Bombastus	Castor oil	110816-3	1
Lamotte/Phönix	Castor oil	1128	1
EuRho	Castor oil	1206039-02	1
EuRho/Phönix	Castor oil	1212008-01	1
Lamotte/gehe	Castor oil	1241892	1
Lamotte/gehe	Castor oil	125162	1
Lamotte/gehe	Castor oil	125473	1
Fagron	Castor oil	126F07-B06-324379	1
Fagron	Castor oil	127170	1
Sanacorp	Castor oil	128224	1
Euro OTC	Castor oil	128224	1
Fagron	Castor oil	13/27-B01-287857	1
Eurho	Castor oil	1304013-01	1
Sanacorp	Castor oil	130260	1
Großhandel	Castor oil	1304013-01	1
EuRho/Geilenkirchen	Castor oil	1304013-02	1
Großhandel	Castor oil	1304013-02	1
Euro OTC	Castor oil	130401302	1
EuRho	Castor oil	1304013-02	1
Caelo	Castor oil	13065401	1
Euro OTC	Castor oil	130401313-02	1
Fagron	Castor oil	13065405	1
EuRho/Noweda	Castor oil	1308019-01	2
Euro OTC	Castor oil	13071301	2
Noweda	Castor oil	1308019-01	1
Phönix	Castor oil	1308019-01	1
Caelo	Castor oil	1308019-02	1
Euro OTC	Castor oil	130801902	1
Anzag 04.03.2014	Castor oil	1308019-02	1
Großhandel	Castor oil	1308019-02	1
Euro OTC	Castor oil	1308019-02	2
Caelo	Castor oil	13134603	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Castor oil	131344604	1
Caelo	Castor oil	13134608	6
Caelo	Castor oil	13174505	2
Caelo	Castor oil	13174507	1
Caelo	Castor oil	1319609	1
Sanacorp	Castor oil	13196308	1
Caelo	Castor oil	13196308	1
Noweda	Castor oil	13196313	1
Fagron	Castor oil	133085	1
Caelo	Castor oil	13357506	1
	Castor oil	13341606	1
Fagron	Castor oil	13G29B01	1
Fagron	Castor oil	13G29-B01-285422	1
Fagron	Castor oil	13I27-B01-287858	1
Caelo	Castor oil	14020101	1
Caelo	Castor oil	13I27-B01-287858	1
Fagron	Castor oil	13L27B01287857	1
Caelo	Castor oil	14020104	2
Sanacorp	Castor oil	14020105	1
Fagron	Castor oil	13I27-B01-287857	1
Caelo	Castor oil	14020105	3
Caelo	Castor oil	14020106	2
Caelo	Castor oil	14020108	11
Fagron	Castor oil	14020108	1
Caelo	Castor oil	14020109	6
AHD/Caealo 15.08.14	Castor oil	14020109	1
Euro OTC	Castor oil	140302401	1
Phönix	Castor oil	1403024-01	1
Euro OTC	Castor oil	1403024-02	8
Caelo	Castor oil	1403024-02	2
Caelo	Castor oil	14058809	1
Euro OTC	Castor oil	140302402	1
Caelo	Castor oil	14058808	1
Euro OTC	Castor oil	140705301	2
AHD 27.11.14 EK 2,53	Castor oil	1407053-01	1
Ph.Eur.8.0	Castor oil	140705301	1
Eur	Castor oil	1407053-01	1
Hedinger / Noweda	Castor oil	1407053-02	3
Sana 10.02.15 EK 2,57	Castor oil	1407053-02	1
Fagron	Castor oil	1407053-02	1
Castor Products Company...	Castor oil	1407053-02	1
Euro OCT	Castor oil	14112006-02	1
Caelo	Castor oil	14115405	1
Euro OTC	Castor oil	141200601	1
Caelo	Castor oil	1412006-01	2
Bombastus	Castor oil	1412006-01	1
Caelo	Castor oil	1412006-02	2
Euro/Noweda	Castor oil	1412006-02	1
Euro OTC	Castor oil	141200602	2
EuRho	Castor oil	141-2006-02	1
Caelo	Castor oil	14181102	1
Caelo	Castor oil	14181103	6
Caelo	Castor oil	14199507	2
Caelo	Castor oil	14199602	1
Caelo	Castor oil	14199604	7
Caelo	Castor oil	14199605	7
Caelo	Castor oil	14199609	8
Caelo	Castor oil	1420106	1
Cealo	Castor oil	14199609	1
Caelo	Castor oil	14304803	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Castor oil	14330601	1
Caelo	Castor oil	14304807	1
Caelo	Castor oil	14304806	1
Caelo	Castor oil	14330601	2
Sanacorp	Castor oil	14330601	1
Caelo	Castor oil	14330604	6
Caelo	Castor oil	14330605	6
Caelo	Castor oil	14330606	10
Fagron	Castor oil	14C17-B01	3
Fagron	Castor oil	14A28-B07-304816	1
Fagron	Castor oil	14C17-B01-293833	2
Fagron	Castor oil	14c17-b01-293832	1
Caelo	Castor oil	15/355/Ri	1
Fagron	Castor oil	14C17-B01-293834	1
Sanacorp	Castor oil	1504006-01	1
Fagron	Castor oil	1504006-01	1
Euro OTC	Castor oil	1504006-01	8
Euro OTC	Castor oil	1504006-02	5
Euro OTC	Castor oil	150400602	1
Caelo	Castor oil	150400602	1
EuRho / 11	Castor oil	1504006-02	2
Unisan	Castor oil	1504006-02	1
Anzag/Henry Lamotte	Castor oil	1504006-02	1
Euro OTC	Castor oil	15072014B	1
Caelo	Castor oil	15067306	1
Phönix Gotha	Castor oil	1509003-01	1
Sanacorp	Castor oil	1509003-01	1
EuRho	Castor oil	1509003-01	1
Euro OTC	Castor oil	150900301	1
Caelo	Castor oil	1509003-01	1
Castor/Phönix	Castor oil	1509003-02	1
Euro OTC	Castor oil	1512004-02	1
Euro OTC	Castor oil	1512024-01	11
EuRho	Castor oil	1512024-01	1
Castor Pruducts Company...	Castor oil	1512024-01	1
Euro OTC	Castor oil	151204-02	1
13.10.2016 Phoenix 3,77EUR	Castor oil	1512024-02	1
Caelo	Castor oil	15144302	5
Castor Pruducts Company...	Castor oil	1512024-02	1
Caelo	Castor oil	15144303	6
sanacorp	Castor oil	15144302	2
Caesar&Loretz/Caesar&Lo...	Castor oil	15144304	1
AllianceHealthcare	Castor oil	15144304	1
Caelo	Castor oil	15144305	2
Caelo	Castor oil	15144307	2
Caelo	Castor oil	15145101	5
Rondell Apotheke	Castor oil	15145101	1
Caelo	Castor oil	15145104	1
Caelo	Castor oil	15145105	1
Caelo	Castor oil	15253404	4
Caelo	Castor oil	15145106	1
Caelo	Castor oil	15253403	1
Fagron	Castor oil	1528-B07-3048	1
Caelo	Castor oil	15392002	2
Caelo	Castor oil	15392003	3
Caelo	Castor oil	15392007	3
Caelo	Castor oil	15399806	1
Euro OTC	Castor oil	154006-02	1
Fagron	Castor oil	15A28B07	1
Fagron	Castor oil	15A28-B07-304818	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Castor oil	15A28-B07-304817	1
Fagron	Castor oil	15A28-B07-304816	2
Fagron	Castor oil	15D02-B02-311714	2
Sanacorp	Castor oil	15D02-Bo2-311714	1
Fagron	Castor oil	15D02-B02-314240	1
Fagron	Castor oil	15J15-B08	4
Fagron	Castor oil	15J15-B08-316099	1
Caelo	Castor oil	16016202	4
Fagron	Castor oil	15J15-B08-318730	1
Caelo	Castor oil	16/081/R	1
Caelo	Castor oil	15J15-B08-318730	1
Anzag	Castor oil	16016207	1
Caelo	Castor oil	16017102	1
Caelo	Castor oil	1604019-01	1
Euro OTC	Castor oil	160401902	1
Caelo	Castor oil	16102801	6
Euro OTC	Castor oil	1604019-02	1
Caelo	Castor oil	16102802	4
Calo/Sanacorp	Castor oil	16102802	2
Caesar & Loretz GmbH	Castor oil	16102802	1
Allianz Healthcare	Castor oil	16102804	1
Henry Lamotte Oils/ San...	Castor oil	16102804	1
Caelo	Castor oil	16102804	2
Euro OTC	Castor oil	161100501	1
Euro OTC	Castor oil	1611005-01	2
Euro OTC	Castor oil	1611005-02	2
Finteler Apotheke	Castor oil	1611005-01	1
Caelo	Castor oil	1611005-01	1
Caelo	Castor oil	16167203	1
Bombastus	Castor oil	161101R/301591	1
Caelo	Castor oil	16200007	4
Caelo	Castor oil	1620007	1
Caesar & Loretz GmbH	Castor oil	16236202	1
Now/Cael030.12.16 EUR2,25	Castor oil	16236202	1
Caelo	Castor oil	16273101	1
Caelo	Castor oil	16236205	1
Caesar & Loretz GmbH	Castor oil	16273105	1
Caelo	Castor oil	16314602	1
Caelo	Castor oil	16314604	2
Caelo	Castor oil	16314603	1
Fagron	Castor oil	16C29-B01	1
Fagron	Castor oil	16C29-B01-320886	1
Fagron	Castor oil	16F07-B06-324379	3
Fagron	Castor oil	16F07-B06	1
Fagron	Castor oil	16H05-B03	2
Fagron	Castor oil	16H05-B03-326864	5
Fagron	Castor oil	16h05-b03-326864	2
Fagron	Castor oil	16J12-B01	1
Hubertus Apotheke	Castor oil	16J12-B01-328159	1
Caelo	Castor oil	17027206	1
Fagron	Castor oil	16J12-B01-328159	2
Caelo	Castor oil	17002280	1
Bombastus	Castor oil	17001711	1
Caelo	Castor oil	17027201	1
Euro OTC	Castor oil	1705035-01	4
Caelo	Castor oil	17107403	3
Caelo	Castor oil	17107501	1
Caelo	Castor oil	17107407	1
Caelo	Castor oil	17107503	2
GEHE	Castor oil	17107506	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Castor oil	17107506	1
Caesar & Loretz GmbH	Castor oil	17113803	1
Caelo	Castor oil	17113802	1
Caelo	Castor oil	17113803	6
Caelo	Castor oil	17189909	1
Caelo	Castor oil	17251901	5
Caelo	Castor oil	17359306	1
Bombastus	Castor oil	18001340	3
Fagron	Castor oil	17H22-B01-341047	1
Caelo	Castor oil	1913A-02009	1
Caelo	Castor oil	18021505	1
Lamotte/Gehe	Castor oil	1914	2
Bombastus	Castor oil	180306R/311684	1
Caelo	Castor oil	20161102-3	1
Caelo	Castor oil	26041403	1
Caelo	Castor oil	24101402	1
Bombastus	Castor oil	268699	1
Gehe 05.03.2014	Castor oil	278660	1
Caelo	Castor oil	28051512	1
Bom	Castor oil	280791	1
	Castor oil	280791	1
Sanacorp	Castor oil	281943	1
Fargon/Sanacorp	Castor oil	281943	4
Bombastu-Werke AG/Gehe	Castor oil	281944	1
Bombastus	Castor oil	283620	7
Gehe	Castor oil	283620	3
Bombastus	Castor oil	284505	6
Caelo	Castor oil	284505	2
Bombastus	Castor oil	285757	3
GEHE 17.07.2014	Castor oil	285757	1
Bombastus	Castor oil	285801	12
Bombastus	Castor oil	287167	12
Bombastus	Castor oil	288126	13
	Castor oil	288126	1
Bombastus/ Anzag	Castor oil	288126	1
Gehe	Castor oil	288126	1
Caelo	Castor oil	288126	1
Bombastus	Castor oil	288638	6
Gehe Kassel	Castor oil	288638	2
Bombastus	Castor oil	288639	16
lamotte Phoenix	Castor oil	288721	1
Gehe Kassel 06.06.2015	Castor oil	290212	1
Bombastus	Castor oil	290213	14
Euro OTC	Castor oil	290714b	1
Bombastus	Castor oil	290714a	1
Hedinger/Sanacorp	Castor oil	290213	1
Anzag	Castor oil	292409	1
Caelo	Castor oil	292409	1
Bombastus	Castor oil	294667	4
Caelo	Castor oil	294667	1
Caelo	Castor oil	294668	1
Fagron	Castor oil	294668	1
Caelo	Castor oil	295412	2
Caelo	Castor oil	296751	3
Bombastus	Castor oil	297429	5
Bombastus	Castor oil	297430	9
Caelo	Castor oil	298698	1
Henry Lamotte	Castor oil	2X1000004554	1
Alliance	Castor oil	300629	1
Caelo	Castor oil	30031504	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Castor oil	300628	1
Caelo	Castor oil	300629	2
Bomastus-Werke/Phönix	Castor oil	300629	1
Bombastus	Castor oil	301587	2
Caelo	Castor oil	301587	1
Bombastus	Castor oil	301591	12
Caelo	Castor oil	301591	3
Fagron	Castor oil	301591	1
Caelo	Castor oil	303857	3
Bombastus	Castor oil	303857	2
Bombastus	Castor oil	303858	10
Euro OTC	Castor oil	303858	1
Bombastus	Castor oil	305229	2
Bombastus	Castor oil	305230	12
phönix	Castor oil	305230	2
Bombastus	Castor oil	307456	5
Bombastus	Castor oil	307457	4
Caelo	Castor oil	307457	3
Bombastus	Castor oil	308825	7
Bombastus	Castor oil	308824	1
Caelo	Castor oil	308825	3
Bombastus	Castor oil	3888	1
Caelo	Castor oil	4614M-02009	1
Bombastus	Castor oil	311684	1
Caelo	Castor oil	4330605	1
Lamotte/Gehe	Castor oil	4198	1
Euro OTC	Castor oil	31032015A	1
Bombastus	Castor oil	4958	2
Defektur	Castor oil	5253	2
Defektur	Castor oil	5052	3
Defektur	Castor oil	5196	1
Defektur	Castor oil	5334	3
Defektur	Castor oil	5408	2
Phönix/Lamotte	Castor oil	5450801	13
Henry Lamotte/ Noweda	Castor oil	5450801	2
Henry Lamotte Oils GmbH	Castor oil	5450801	1
Lamotte/L3	Castor oil	5450801	1
Henry Lamotte Oils, Gehe	Castor oil	5450801	1
Henry Lamotte/Gehe	Castor oil	5450801	2
Fagron	Castor oil	5450801	1
Lamotte/Alliance	Castor oil	5450801	1
Rez	Castor oil	5519	2
Euro OTC	Castor oil	5581	1
Defektur	Castor oil	5775	1
Henry Lamotte Oils/GEHE	Castor oil	5899201	1
Noweda/Henry Lamotte	Castor oil	5899201	2
Henry Lamotte/Sanacorp	Castor oil	5899201	2
Henry la Motte	Castor oil	5899201	1
Henning/ Phoenix	Castor oil	5899201	2
lamotte/Phönix	Castor oil	5899201	1
Caelo	Castor oil	5899201	1
Defektur Abfüllung	Castor oil	5945	1
Defektur Abfüllung	Castor oil	6001	2
Euro OTC	Castor oil	60082103	1
Bombastus	Castor oil	60082103	2
Defektur Abfüllung	Castor oil	6081	1
Defektur Abfüllung	Castor oil	6260	1
Defektur Abfüllung	Castor oil	6189	1
Defektur Abfüllung	Castor oil	6333	1
Defektur	Castor oil	6442	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Defektur Abfüllung	Castor oil	6414	1
Bombastus	Castor oil	91696	1
Euro OTC	Castor oil	ChB120115	1
Caelo	Castor oil	Fertigarzneimittel	1
ACM	Castor oil	KNC0899J15	1
Ichthyol	Castor oil		1
Fagron	Castor oil		1
Caelo	Castor oil	1206039-01	1
Eurho / Sanacorp	Castor oil	121008-02	1
Caelo	Castor oil	11359507	1
Caelo	Castor oil	1/AKD1	1
Euro OTC	Castor oil	126279	1
Fagron	Castor oil	12D04-N09	2
Caelo	Castor oil	1304013-01	1
Cealo	Castor oil	13096304	1
Bergchemie	Castor oil	13096304	1
Fagron	Castor oil	13D29-N01	1
Euro OTC	Castor oil	1513A03553	2
Euro OTC	Castor oil	3313U-03553	1
Euro OTC	Castor oil	3513I-03553	1
Bombastus	Castor oil	268877	1
Bombastus	Castor oil	271753	2
Bombastus	Castor oil	273389	1
Caelo	Castor oil	272591	1
Noweda	Castor oil	273390	1
Caelo	Castor oil	269565	1
Phoenix 21.03.13	Castor oil	273881	2
Phönix 07.12.12	Castor oil	273390	1
Bombastus	Castor oil	273881	2
Henry Lamotte Oils/Gehe	Castor oil	123384	1
Caelo	Castor oil	11359508	1
Euro OTC	Castor oil	1202041-01	1
Geo-Chem/sanacorp	Castor oil	1206039-01	1
Euro OTC	Castor oil	1206039-02	1
Caelo	Castor oil	12071107	1
Geo-Chem/sanacorp	Castor oil	1212008-01	1
Sanacorp	Castor oil	12149306	1
Noweda	Castor oil	1304013-01	1
Caelo	Castor oil	15012016	1
gehe, 18.01.2013	Castor oil	272591	1
Gehe	Castor oil	275937	1
Gehe	Castor oil	275938	2
Caelo	Castor oil	275938	3
Bombastus	Castor oil	278660	3
Lamotte	Castor oil	5188401	1
Noweda / Henry Lamotte	Castor oil	5280301	1
HenryLamotte/Phoenix	Castor oil	5280301	1
Otto Fischar Gmbh	Castor oil	5334001	1
Henry Lamotte	Castor oil	5334001	1
Lamotte	Castor oil	5334001	2
lamotte	Castor oil	5450801	1
lamotte/anzag	Castor oil	5334001	1
Henry Lamotte Oils	Castor oil	5334001	1
Ph.Eur.	Castor oil	6861795	1

- 89 775 spectra from 501 *Apo-Ident* customers from a total of 26 700 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are docu-

mented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Castor oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Castor oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	285	0	35 259
Type B	0	795	0	93 079
Type C	0	1574	3	89 662

The substance/substance group *Castor oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.8947 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 99.2453 %)
Type C	100.0000 % (> 99.6586 %)	99.8098 % (> 99.6195 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30741	31343	2.03	47.98
31343	31343	0.00	50.19
31739	31739	0.00	47.18
31751	31739	1.69	47.56
31824	31739	0.79	46.70
31845	31845	0.00	45.70
32943	31739	4.22	48.54

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Cedar leaf oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31106-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Cedar leaf oil; Oleum cedri e fol. aether.

Special notes

When selecting the *Cedar leaf oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cedar leaf oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cedar leaf oil	13337901	31106	40	1404622
Caelo	Cedar leaf oil	13337901	31200	40	20131011
Caelo	Cedar leaf oil	13337902	31421	40	20131011

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Cedar leaf oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 424 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 72 spectra from 3 reference samples of the substance/substance group *Cedar leaf oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cedar leaf oil	13337902	31421 [†]	20
Caelo	Cedar leaf oil	15305207	33607	40
Caelo	Cedar leaf oil	15305206	not required	12

- 93 802 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 0 spectra from 0 *Apo-Ident* customers from 0 batches from the substance/substance group *Cedar leaf oil*.
- These include spectra of independent samples from 0 batches, of which no spectra have been used for the generation of the database.
- 91 352 spectra from 501 *Apo-Ident* customers from a total of 27 009 batches from a further 334 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cedar leaf oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cedar leaf oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	72	0	93 802
Type C	0	0	0	91 352

The substance/substance group *Cedar leaf oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 91.6667 %)
Type C	100.0000 % (> 99.6586 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31106	31106	0.00	87.67
31200	31106	1.36	88.55
31421	31106	2.07	88.33

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Cedar wood oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30938-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Cedar wood oil; Oleum cedri ligni aether.

Special notes

When selecting the *Cedar wood oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cedar wood oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cedar wood oil	12180708	30938	50	1404506
Caelo	Cedar wood oil	12180710	31267	40	20120628
Caelo	Cedar wood oil	14230204	32455	40	1610271

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 130 spectra of 3 reference samples from the substance/substance group *Cedar wood oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 414 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 80 spectra from 2 reference samples of the substance/substance group *Cedar wood oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cedar wood oil	14230211	32748	40
Caelo	Cedar wood oil	16028815	33785	40

- 93 794 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 9 spectra from 9 *Apo-Ident* customers from 7 batches from the substance/substance group *Cedar wood oil*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Cedar wood oil	12180709	1
Caelo	Cedar wood oil	12180710	1
Caelo	Cedar wood oil	14230211	3
Caelo	Cedar wood oil	16028805	1
Caelo	Cedar wood oil	16028815	1
Phönix	Cedar wood oil	308062	1
Caelo	Cedar wood oil	12180703	1

- 91 343 spectra from 501 *Apo-Ident* customers from a total of 27 002 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cedar wood oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cedar wood oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	130	0	35 414
Type B	0	80	0	93 794
Type C	0	8	1	91 343

The substance/substance group *Cedar wood oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.3846 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 92.5000 %)
Type C	100.0000 % (> 99.6595 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30938	30938	0.00	25.01
31267	30938	1.90	25.36
32455	32455	0.00	30.64

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Cinnamon bark oil, artificial**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30605-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Cinnamon bark oil, artificial; Oleum Cinnamomi artificiale

Special notes

When selecting the *Cinnamon bark oil, artificial* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cinnamon bark oil, artificial*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cinnamon bark oi...	11275903	30605	40	1402333
Caelo	Cinnamon bark oi...	14377301	31617	40	20150107

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Cinnamon bark oil, artificial*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Cinnamon bark oil, artificial*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cinnamon bark oil, artificial	14377301	31617 [†]	20
Caelo	Cinnamon bark oil, artificial	16026505	33609	40
Caelo	Cinnamon bark oil, artificial	16026505	33610	40
Caelo	Cinnamon bark oil, artificial	180217	34057	40

- 93 734 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 6 spectra from 6 *Apo-Ident* customers from 6 batches from the substance/substance group *Cinnamon bark oil, artificial*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Cinnamon bark oil, artificial	2031501	1
Caelo	Cinnamon bark oil, artificial	14377304	1
Caelo	Cinnamon bark oil, artificial	14377307	1
Caelo	Cinnamon bark oil, artificial	11275903	1
Caelo	Cinnamon bark oil, artificial	10112014	1
Caelo	Cinnamon bark oil, artificial	12019604	1

- 91 346 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cinnamon bark oil, artificial* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cinnamon bark oil, artificial* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	140	0	93 734
Type C	0	5	1	91 346

The substance/substance group *Cinnamon bark oil, artificial* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6600 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30605	30605	0.00	21.59
31617	30605	3.67	20.89

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Cinnamon bark oil, Ceylon
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30251-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Cinnamon bark oil, Ceylon; Oleum cinnamomi ceylanici

Special notes

When selecting the *Cinnamon bark oil, Ceylon* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cinnamon bark oil, Ceylon*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cinnamon bark oi...	11352506	30604	40	20111219
Caelo	Cinnamon bark oi...	14092802	31560	40	AR-15-FG-001171-01
Caelo	Cinnamon bark oi...	163243	33091	40	20161129*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Cinnamon bark oil, Ceylon*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Cinnamon bark oil, Ceylon*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cinnamon bark oil, Ceylon	10185201	30251	40
Caelo	Cinnamon bark oil, Ceylon	14092802	31560 [†]	20
Caelo	Cinnamon bark oil, Ceylon	15278210	33480	40
Caelo	Cinnamon bark oil, Ceylon	15278218	33783	40

- 93 734 spectra from a total of 2163 batches from further 436 substances. These spectra were

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 42 spectra from 25 *Apo-Ident* customers from 21 batches from the substance/substance group *Cinnamon bark oil, Ceylon*.
- These include spectra of independent samples from 19 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Ölbaum / Anzag	Cinnamon bark oil, Ceylon	10220032803	1
Apotheke in Steinkuhl	Cinnamon bark oil, Ceylon	110504/3	1
Caelo	Cinnamon bark oil, Ceylon	11352506	1
Caelo	Cinnamon bark oil, Ceylon	13044802	5
Caelo	Cinnamon bark oil, Ceylon	12181408	1
Taoasis	Cinnamon bark oil, Ceylon	130111-110762	2
/Noweda	Cinnamon bark oil, Ceylon	13044802	4
Caelo	Cinnamon bark oil, Ceylon	13044810	5
Caelo	Cinnamon bark oil, Ceylon	15278206	3
Caelo	Cinnamon bark oil, Ceylon	11352501	1
Taoasis GmbH	Cinnamon bark oil, Ceylon	130111-110345	1
Dr.Schlegels Resana/NOW	Cinnamon bark oil, Ceylon	1205-3566	1
Caelo	Cinnamon bark oil, Ceylon	14071203	1
Caelo	Cinnamon bark oil, Ceylon	14092802	3
Caelo	Cinnamon bark oil, Ceylon	13044808	1
Caelo	Cinnamon bark oil, Ceylon	14092805	5
Dr.Schlegels Resana/NOW	Cinnamon bark oil, Ceylon	1407-4584	1
Caelo	Cinnamon bark oil, Ceylon	13044809	1
Caelo	Cinnamon bark oil, Ceylon	14092806	1
Caelo	Cinnamon bark oil, Ceylon	15151106	1
Caelo	Cinnamon bark oil, Ceylon	16324307	1
Caelo	Cinnamon bark oil, Ceylon	12181403	1

- 91 310 spectra from 501 *Apo-Ident* customers from a total of 26 988 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cinnamon bark oil, Ceylon* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cinnamon bark oil, Ceylon* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	140	0	93 734
Type C	0	38	4	91 310

The substance/substance group *Cinnamon bark oil, Ceylon* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6588 %)	90.4762 % (> 83.3333 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30604	31560	5.72	32.40
31560	31560	0.00	33.87
33091	33091	0.00	27.44

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Citronella oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30208-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Citronella oil; Oleum cymbopogon winterianus

Special notes

When selecting the *Citronella oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Citronella oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Citronella oil	12359903	31251	40	1405526
Caelo	Citronella oil	13287808	31408	40	20130827
Caelo	Citronella oil	14229803	31631	40	20140825

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Citronella oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 181 spectra from 5 reference samples of the substance/substance group *Citronella oil*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Citronella oil	09323007	30208	41
Caelo	Citronella oil	91668379	30250	40
Caelo	Citronella oil	14229803	31631 [†]	20
Caelo	Citronella oil	162791	32692	40
Caelo	Citronella oil	172831	33675	40

- 93 693 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 47 spectra from 33 *Apo-Ident* customers from 38 batches from the substance/substance group *Citronella oil*.
- These include spectra of independent samples from 35 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Citronella oil	12359903	2
Caelo	Citronella oil	12359907	1
Caelo	Citronella oil	13287802	1
Caelo	Citronella oil	13287808	1
Caelo	Citronella oil	13287809	1
Caelo	Citronella oil	14229803	1
Caelo	Citronella oil	14229805	1
Caelo	Citronella oil	15222302	1
Caelo	Citronella oil	15257304	2
Caelo	Citronella oil	286270	1
Bombastus	Citronella oil	286270	2
Phönix Gotha	Citronella oil	288826	1
Bergland	Citronella oil	C400822	1
Caelo	Citronella oil	15222303	1
Caelo	Citronella oil	15257309	1
Caelo	Citronella oil	15257308	1
Bombastus	Citronella oil	297193	1
Caelo	Citronella oil	15257305	1
Caelo	Citronella oil	16279105	1
Bombastus	Citronella oil	297194	1
Bombastus	Citronella oil	302910	1
Bombastus	Citronella oil	303977	1
Bombastus	Citronella oil	302908	1
Caelo	Citronella oil	303977	2
Bombastus	Citronella oil	309032	1
Caelo	Citronella oil	10193312	1
Caelo	Citronella oil	11074203	1
Caelo	Citronella oil	11074201	1
Bombastus	Citronella oil	1213A-01974	1
Caelo	Citronella oil	11101301	1
Caesar/Lorenz	Citronella oil	12224603	1
Euro OTC	Citronella oil	12224603	1
Caelo	Citronella oil	12359902	1
Kaders	Citronella oil	131097	1
Caelo	Citronella oil	13287804	2
Caelo	Citronella oil	12359909	1
Phönix Gotha	Citronella oil	260477	1
Bombastus	Citronella oil	269807	1
Bombastus	Citronella oil	276241	1
Bombastus	Citronella oil	280657	1
Caelo	Citronella oil	276241	1
Taoasis	Citronella oil	104719-125577	1

- 91 305 spectra from 501 *Apo-Ident* customers from a total of 26 972 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are docu-

mented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Citronella oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Citronella oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	154	27	93 693
Type C	0	45	2	91 305

The substance/substance group *Citronella oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	85.0829 % (> 83.4254 %)
Type C	100.0000 % (> 99.6587 %)	95.7447 % (> 89.3617 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31251	31251	0.00	42.25
31408	31251	6.98	42.37
31631	31251	2.47	43.28

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Citrus oils
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31017-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Citrus oils; Bitter orange oil; Citric oil; Grapefruit oil; Limonis aetheroleum; Oleum aurantii dulcis; Oleum aurantii pericarpium; Oleum citri; Oleum citri paradisi artificiale; Sweet orange oil

Special notes

When selecting the *Citrus oils* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Citrus oils*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Grapefruit oil	12367403	31017	40	1404621
Fagron	Citric oil	13E06-N10	31046	40	20130618

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sweet orange oil	12043412	31164	40	1405105
Caelo	Sweet orange oil	13378901	31206	40	1405107
Fagron	Citric oil	12G13-B06-287454	31261	40	1407356
Caelo	Sweet orange oil	13024204	31273	40	20130124
Caelo	Sweet orange oil	13378906	31295	40	1407188
Caelo	Sweet orange oil	13401605	31297	38	20131220
Caelo	Sweet orange oil	13024204	31411	40	20130124
Caelo	Grapefruit oil	13252306	31417	40	20130802
Caelo	Citric oil	13259522	31432	40	20130816
Euro OTC	Citric oil	1202024-01	31524	40	20120312
Caelo	Bitter orange oil	13303803	31861	40	1606482
Caelo	Bitter orange oil	15397901	32181	40	20151125
Caelo	Bitter orange oil	170369	33062	30	20170224*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 588 spectra of 15 reference samples from the substance/substance group *Citrus oils*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 14 different batches.
- 34 956 spectra from a total of 826 batches from further 208 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 560 spectra from 16 reference samples of the substance/substance group *Citrus oils*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Sweet orange oil	09319501	30171	40
Caelo	Citric oil	90965219	30248	20
Caelo	Grapefruit oil	13252306	31417 [†]	20
Caelo	Sweet orange oil	13378906	31418	60
Caelo	Citric oil	13259522	31432 [†]	20
Euro OTC	Citric oil	1202024-01	31524 [†]	20
Caelo	Bitter orange oil	13303803	31861 [†]	20
Caelo	Bitter orange oil	15397901	32351	40
Caelo	Citric oil	162752	32687	40
Caelo	Citric oil	171068	33305	40
Caelo	Sweet orange oil	171596	33307	40
Caelo	Sweet orange oil	16015501	33467	40
Caelo	Grapefruit oil	15429104	33630	40
Caelo	Grapefruit oil	172307	34054	40
Caelo	Bitter orange oil	180889	34078	40
Caelo	Grapefruit oil	182318	34237	40

- 93314 spectra from a total of 2151 batches from further 433 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 344 spectra from 146 *Apo-Ident* customers from 165 batches from the substance/substance group *Citrus oils*.
- These include spectra of independent samples from 155 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Taoasis	Citric oil	4LI14M-118369	1
Taoasis	Citric oil	122527-122620	1
Fagron	Citric oil	11C28-N060363/440	2
Caelo	Citric oil	12360012	8
Caelo	Citric oil	12360002	6
Caelo	Citric oil	12367408	1
Fagron	Citric oil	12G13-B06-287454	1
Caelo	Citric oil	13259508	8
Caelo	Citric oil	13259504	8
Caelo	Citric oil	13259523	6
Fagron	Citric oil	13B01N02	2
Fagron	Citric oil	13K14-B03-294706	5
Fagron	Citric oil	13K14-B03-289819	1
Fagron	Citric oil	13E06-N10	6
Caelo	Citric oil	14121807	5
Taoasis	Citric oil	1401270-121612	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Citric oil	13K14-B03-787033	1
Taoasis	Citric oil	140116-115338	1
Caelo	Citric oil	14121808	2
Caelo	Citric oil	14121815	5
Caelo	Citric oil	14121813	1
Caelo	Citric oil	14365201	3
Caelo	Citric oil	14365202	4
Caelo	Citric oil	14365212	14
Fagron	Citric oil	14H06B02308273	1
Caelo	Citric oil	14H06-B02	1
Fagron	Citric oil	14H06-B02-299949	2
Fagro/Gehe	Citric oil	14H06-B02-299949	1
Caelo	Citric oil	15182702	1
Fagron	Citric oil	14H06-B02-308273	2
Caelo	Citric oil	15121403	1
Caelo	Citric oil	15182703	11
Caelo	Citric oil	288778	1
Bombastus	Citric oil	280999	4
Bombastus	Citric oil	289673	1
Bombastus	Citric oil	294337	2
Sanacorp	Citric oil	5/AKD1	1
Fagron	Citric oil	733117	1
Bombastus	Citric oil	294951	1
Bombastus	Citric oil	284045	1
Sanacorp	Citric oil		1
Caelo	Citric oil	202170C16028702	1
Taoasis	Citric oil	3J16M-125453	1
Taoasis	Citric oil	5LI12M-109229	1
Caelo	Citric oil	9322802	1
Caelo	Citric oil	9322807	1
Noweda	Citric oil	10375003	2
Caelo	Citric oil	9322815	1
Caelo	Citric oil	10375003	2
Caelo	Citric oil	10375012	1
Caelo	Citric oil	10375014	2
Caelo	Citric oil	11231807	2
Fagron	Citric oil	11D28-N09	1
Caelo	Citric oil	1135408	1
Fagron	Citric oil	11B04-N02	2
Fagron	Citric oil	11K28-N04	2
Fagron	Citric oil	11I05-N13	1
DEFEKTUR	Citric oil	1202024-01	1
Fagron	Citric oil	11K28-N04	1
Caelo	Citric oil	12040502	4
Caelo	Citric oil	12040504	9
Sanacorp	Citric oil	12040504	1
Anzag 15.11.12 ek:17,18	Citric oil	12040507	1
Caelo	Citric oil	12040507	1
Taoasis	Citric oil	122650-123634	1
Caelo	Citric oil	123600002	1
Caelo	Citric oil	123600001	2
Caelo	Citric oil	123600001	1
Anzag 17.06.2013	Citric oil	123600001	1
Noweda	Citric oil	123600002	1
Caelo	Citric oil	123600010	2
Caelo	Citric oil	12360010c	1
Caelo	Citric oil	12360011	11
Caelo	Citric oil	123600002	1
Caelo	Citric oil	12360013	1
Caelo	Citric oil	12360012x2	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Citric oil	12A23-N05	7
Fagron	Citric oil	12A23-N05	1
Fagron	Citric oil	12G04-N21	2
Caelo	Citric oil	1325904	1
Caelo	Citric oil	132259504	1
Caelo	Citric oil	1359512	1
Gehe 29.04.2013	Citric oil	13B01N/02	1
Fagron	Citric oil	13B01-N02	6
Fagron	Citric oil	13b01-n02	1
Caelo	Citric oil	13E06-N100875/900	1
Caelo	Citric oil	15051314	1
Caelo	Citric oil	15182711	4
Sanacorp	Citric oil	15182711	1
Fagron	Citric oil	15F29-B04-309901	1
Caelo	Citric oil	15182718	2
Fagron	Citric oil	15F29-B04-322985	2
Fagron	Citric oil	15F29-B04	2
Fagron	Citric oil	15F29-B04-327600	1
Caelo	Citric oil	16028702	4
Caelo	Citric oil	16028715	1
Caelo	Citric oil	16028709	1
Caelo	Citric oil	16028719	6
Caelo	Citric oil	16028712	1
Caelo	Citric oil	16275203	2
Caelo	Citric oil	16275202	1
Caelo	Citric oil	16275216	2
Taoasis	Citric oil	1786-125280	1
Caelo	Citric oil	20111307	1
Caelo	Citric oil	2515I-01973	1
Fagron	Citric oil	268961	1
Caelo	Citric oil	28051308	1
Bombastus	Citric oil	266679	1
Bombastus	Citric oil	297734	1
Bombastus	Citric oil	299753	1
Taoasis	Citric oil	2LI08M-73396	1
Caelo	Citric oil	3710A-01973	1
Caelo	Citric oil	31121206	1
Bombastus	Citric oil	300294	1
Taoasis	Citric oil	3805201	1
Fagron	Citric oil	602796	1
Fagron	Citric oil	645920	4
Fagron	Citric oil	709110	1
Caelo	Grapefruit oil	13252303	1
Caelo	Grapefruit oil	13252306	1
Caelo	Grapefruit oil	14107304	1
Caelo	Grapefruit oil	15171201	1
Caelo	Grapefruit oil	15429101	2
Caelo	Grapefruit oil	15429105	2
Taoasis	Grapefruit oil	31564-120540	1
Taoasis	Grapefruit oil	44040A-124468	1
Taoasis	Grapefruit oil	333555-111392	1
Apotheker-Bauer	Grapefruit oil		1
Fagron	Citric oil	13K14-B03	2
Bombastus	Citric oil	285568	1
Caelo	Sweet orange oil	15045101	4
Caelo	Sweet orange oil	13024204	4
Caelo	Grapefruit oil	12367403	1
Fagron	Citric oil	12g13-b06-287454	1
VDL;4,32EUR;25.04.15	Citric oil	14121815	1
Fagron	Citric oil	14H06-B02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Citric oil	13259504	1
Bombastus	Citric oil	256839	1
Bombastus	Citric oil	269281	1
Bombastus	Citric oil	284351	1
Caelo	Bitter orange oil	15397903	1
Nachprüfung	Sweet orange oil	12043403	5
Caelo	Sweet orange oil	12043406	1
Caelo	Sweet orange oil	12043409	1
Caelo	Sweet orange oil	13378906	1
Caelo	Sweet orange oil	13389406	3
Taoasis	Sweet orange oil	140129-115917	1
Caelo	Sweet orange oil	14133505	1
Caelo	Sweet orange oil	14150105	1
Caelo	Sweet orange oil	14348803	3
Caelo	Sweet orange oil	14348809	1
Caelo	Sweet orange oil	15051701	1
Caelo	Sweet orange oil	15182904	1
Bombastus	Sweet orange oil	285361	1
Taoasis GmbH	Sweet orange oil	128129	1
Caelo	Sweet orange oil	15285	1
Caelo	Sweet orange oil	15316601	2
Caelo	Sweet orange oil	15402404	1
Caelo	Sweet orange oil	7051301	1
Primavera	Sweet orange oil	530914	1
Phoenix	Sweet orange oil	16015502	1
Caelo	Sweet orange oil	8121204	1
Caelo	Sweet orange oil	16015501	1
Caelo	Sweet orange oil	10375113	1
Fiebig 09.11.12 ek:6,92	Sweet orange oil	11031802	1
Anzag	Sweet orange oil	11310411	1
Großhandel	Sweet orange oil	12043403	2
Caelo	Sweet orange oil	12043410	3
Caelo	Sweet orange oil	12322603	1
Caelo	Sweet orange oil	12322607	1
Caelo	Sweet orange oil	1302404	1
Caelo	Sweet orange oil	13024202	1
Caelo	Sweet orange oil	13303803	1
Caelo	Sweet orange oil	13389401	1
Caelo	Sweet orange oil	15051313	1
Caelo	Sweet orange oil	13389407	1
Phönix 16.05.2013	Sweet orange oil	271897	1
Bombastus	Sweet orange oil	271897	1

- 91 008 spectra from 501 *Apo-Ident* customers from a total of 26 844 batches from a further 330 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Citrus oils* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Citrus oils* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table

breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	588	0	34 956
Type B	0	559	1	93 314
Type C	0	322	22	91 008

The substance/substance group *Citrus oils* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.9796 %)
Type B	100.0000 % (> 99.9761 %)	99.8214 % (> 99.2857 %)
Type C	100.0000 % (> 99.6586 %)	93.6047 % (> 92.7326 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31017	31017	0.00	43.05
31046	31261	9.94	36.41
31164	31164	0.00	43.20
31206	31206	0.00	46.07
31261	31261	0.00	31.35
31273	31206	2.58	49.51
31295	31295	0.00	46.84
31297	31295	2.52	50.14
31411	31295	1.76	48.67
31417	31164	2.35	46.28
31432	31261	2.07	31.01
31524	31261	5.14	31.05
31861	31861	0.00	30.55
32181	31861	8.53	39.84
33062	33062	0.00	39.25

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Clotrimazole 10% Cordes® RK
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30213-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Clotrimazole 10% Cordes® RK

Special notes

When selecting the *Clotrimazole 10% Cordes® RK* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Clotrimazole 10% Cordes® RK*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Ichthyol	Clotrimazole 10%...	11D044	30729	40	20111213
Ichthyol	Clotrimazole 10%...	12D016	31101	40	1507525
Ichthyol	Clotrimazole 10%...	14D009	31806	40	20140606

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Clotrimazole 10% Cordes® RK*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 150 spectra from 4 reference samples of the substance/substance group *Clotrimazole 10% Cordes® RK*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Ichthyol	Clotrimazole 10% Cordes® RK	09D036	30213	50
Ichthyol	Clotrimazole 10% Cordes® RK	14D009	31806 [†]	20
Ichthyol	Clotrimazole 10% Cordes® RK	14D009	32950	40
Ichthyol	Clotrimazole 10% Cordes® RK	16D003	33498	40

- 93 724 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 77 spectra from 18 *Apo-Ident* customers from 9 batches from the substance/substance group *Clotrimazole 10% Cordes® RK*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Ichthyol-Gesellschaft/N...	Clotrimazole 10% Cordes® RK	110044	1
Ichthyol/Alliance Healt...	Clotrimazole 10% Cordes® RK	120016	1
Ichthyol/Noweda	Clotrimazole 10% Cordes® RK	12D016	3
Ichthyol/Jenne	Clotrimazole 10% Cordes® RK	12D016	11
Ichthyol	Clotrimazole 10% Cordes® RK	12D016	8
Ichthyol/AHD	Clotrimazole 10% Cordes® RK	12D016	1
Ichthyol-Gesellschaft/ N...	Clotrimazole 10% Cordes® RK	14D009	2
AHD/ichthyol	Clotrimazole 10% Cordes® RK	14D0089	1
Noweda	Clotrimazole 10% Cordes® RK	14D009	1
ICHTHYOL-GESELLSCHAFT C...	Clotrimazole 10% Cordes® RK	14D009	1
WE 29.11.16 Ichthyol/Sa...	Clotrimazole 10% Cordes® RK	14D009	1
Ichthyol4D/ Noweda	Clotrimazole 10% Cordes® RK	14D009	1
Ichthyol/Noweda	Clotrimazole 10% Cordes® RK	14D009	1
Ichthyol Gesellschaft C...	Clotrimazole 10% Cordes® RK	11D044	1
Cordes	Clotrimazole 10% Cordes® RK	11D044	1
Ichthyol	Clotrimazole 10% Cordes® RK	11D044	25
Ichthyol	Clotrimazole 10% Cordes® RK	11D044	1
Ichthyol Gesellschaft C...	Clotrimazole 10% Cordes® RK	11D044	1
Ichthyol/Jenne	Clotrimazole 10% Cordes® RK	11D044	1
Ichthyol-Gesellschaft/G...	Clotrimazole 10% Cordes® RK	12D016	2
Ichthyol/ Noweda	Clotrimazole 10% Cordes® RK	12D016	1
Ichthyol-Gesellschaft / K	Clotrimazole 10% Cordes® RK	12D016	1
Ichthyol Noweda	Clotrimazole 10% Cordes® RK	12D016	1
Ichthyol	Clotrimazole 10% Cordes® RK	9D036	9

- 91 275 spectra from 501 *Apo-Ident* customers from a total of 27 000 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Clotrimazole 10% Cordes® RK* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Clotrimazole 10% Cordes® RK* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	150	0	93 724
Type C	0	77	0	91 275

The substance/substance group *Clotrimazole 10% Cordes® RK* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.0000 %)
Type C	100.0000 % (> 99.6587 %)	100.0000 % (> 92.2078 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30729	31101	4.94	91.93
31101	31101	0.00	89.52
31806	31101	18.55	103.71

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Clove oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30858-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Clove oil; Caryophylli floris aetheroleum; Oleum caryophylli

Special notes

When selecting the *Clove oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Clove oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Clove oil	12165905	30858	40	20120627
Caelo	Clove oil	13334902	31299	40	20131009
Caelo	Clove oil	13334902	31409	40	20131009
Euro OTC	Clove oil	1210057-01	31526	40	AR-14-FG-013260-01
Euro OTC	Clove oil	1505022	31707	40	0000010252

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Clove oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 344 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 200 spectra from 6 reference samples of the substance/substance group *Clove oil*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Clove oil	1210057-01	31526 [†]	20
Euro OTC	Clove oil	1505022	31707 [†]	20
Caelo	Clove oil	15402108	33299	40
Caelo	Clove oil	171965	33683	40
Euro OTC	Clove oil	1710036	33847	40
Euro OTC	Clove oil	1802017	33999	40

- 93 674 spectra from a total of 2161 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 434 spectra from 140 *Apo-Ident* customers from 130 batches from the substance/substance group *Clove oil*.
- These include spectra of independent samples from 126 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Clove oil	14243105	1
Caelo	Clove oil	14243104	17
Caelo	Clove oil	15119303	2
Caelo	Clove oil	11231706	8
Caelo	Clove oil	13334905	7
Caelo	Clove oil	12380601	10
Fischar/Gehe	Clove oil	13334905	1
Caelo	Clove oil	13334902	14
Caelo	Clove oil	14133602	7
Caelo	Clove oil	15316103	10
Caelo	Clove oil	10272002	3
Caelo	Clove oil	11231702	8
Caelo	Clove oil	11231704	4
Caelo	Clove oil	12165902	6
phönix 30.07.2013	Clove oil	12165905	1
Caelo	Clove oil	13106501	14
Caelo	Clove oil	12165905	14
Caelo	Clove oil	13106502	4
Caelo	Clove oil	12165905	1
Caelo	Clove oil	13182301	21
Caelo	Clove oil	15119304	19
Caelo	Clove oil	15316102	35
Caelo	Clove oil	15316101	2
Fagron	Clove oil	204211969	1
Caelo	Clove oil	2081402	1
Fagron	Clove oil	208211969	1
Caelo	Clove oil	2091306	1
Caelo	Clove oil	4061404	1
Caelo	Clove oil	9308768	1
Caelo	Clove oil	10272012	3
Caelo	Clove oil	6071501	1
Caelo	Clove oil	10272014	2
Euro OTC	Clove oil	1108015-01	2
Anzag 03.07.2012	Clove oil	1108015-01	1
Noweda	Clove oil	11231702	2
Phoenix, 31.08.2012	Clove oil	11231702	1
Caelo	Clove oil	11231703	5
Anzag am 03.07.2012	Clove oil	11231703	1
11231703	Clove oil	11231703	1
Phönix 21.12.12	Clove oil	11K17-N07	1
Caelo	Clove oil	11231709	3
Fagron	Clove oil	11L21-N07	2
Fagron	Clove oil	11K17-N07	1
phönix 22.01.2013	Clove oil	11L21-N07	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Clove oil	1206012-01	1
Phönix, 21.12.2012	Clove oil	11L21N07	1
Caelo	Clove oil	12011502	1
Euro OTC	Clove oil	1206013-01	8
Euro OTC	Clove oil	1210057-01	1
EuRho	Clove oil	1210057-01	2
Euro OTC	Clove oil	121005701	1
Caelo	Clove oil	121655905	1
Euro/Gehe	Clove oil	1210057-01	1
Caesar/Noweda	Clove oil	12165902	1
Noweda	Clove oil	12165905	3
Sanacorp	Clove oil	12165905	1
Phönix	Clove oil	12165905	1
Caelo	Clove oil	12272004-1-0	1
phönix 14.08.2013	Clove oil	12C27-N04	1
Caelo	Clove oil	123806012x	1
Fagron	Clove oil	12C27-N04	1
Fagron	Clove oil	12F01-N01	2
gehe, 10.01.2013	Clove oil	130110	1
Taoasis	Clove oil	130301-110695	2
Calo	Clove oil	13106501	1
Gehe	Clove oil	13106501	1
phönix 30.07.2013	Clove oil	13182301	1
Noweda	Clove oil	13334902	1
Phonix	Clove oil	13334905	1
Caesar/Noweda	Clove oil	13334905	1
Fagron	Clove oil	13L03-B01	2
Fagron	Clove oil	13L03-B01294950	1
Caelo	Clove oil	14133603	4
Caelo	Clove oil	14133608	8
Caelo	Clove oil	14243101	5
Caelo	Clove oil	14243102	1
Caelo	Clove oil	14243103	11
Euro OTC	Clove oil	1505022-02	1
Caelo	Clove oil	15051508	1
Caelo	Clove oil	15101401	1
Caelo	Clove oil	15101402	1
Fagron	Clove oil	1512211969	1
Caelo	Clove oil	1519304	1
Caelo	Clove oil	15316012	1
Phoenix	Clove oil	15316102	1
WE:Sanacorp 13.08.2016...	Clove oil	15316103	1
Caelo	Clove oil	15402101	4
Euro OTC	Clove oil	15402101	1
Caesar & Loretz	Clove oil	15402102	1
Caelo	Clove oil	15402102	1
Caelo	Clove oil	15402108	12
GEHE	Clove oil	15402108	1
Caelo	Clove oil	15402109	2
Sanacorp EK:17,64 EU...	Clove oil	15402109	1
Bombastus	Clove oil	1612020C/296428	1
Caelo	Clove oil	16131401	4
Caelo	Clove oil	16121305	1
Caelo	Clove oil	16131412	8
Caelo	Clove oil	16131413	3
AHD	Clove oil	16131413	1
Caelo	Clove oil	16131416	1
Vögele	Clove oil	2017107539	1
Euro OTC	Clove oil	171003601	1
Caelo	Clove oil	17051402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Clove oil		1
Caelo	Clove oil	20021308	1
Caelo	Clove oil	18071302	1
Vögele	Clove oil	2018100124	1
Caelo	Clove oil	18121411	1
Vögele	Clove oil	2018102519	1
Caelo	Clove oil	265130	1
Caelo	Clove oil	23041404	1
Caelo	Clove oil	23051405	1
Bombastus	Clove oil	265673	1
Bombastus	Clove oil	267699	1
Caelo	Clove oil	23081402	1
Bombastus	Clove oil	272735	1
Caelo	Clove oil	272735	2
Caelo	Clove oil	267700	1
Bombastus	Clove oil	273027	1
Großhandel	Clove oil	273462	2
Bombastus	Clove oil	274495	1
Bombastus	Clove oil	279262	2
Caelo	Clove oil	28051510	1
Euro OTC	Clove oil	274495	1
Bombastus	Clove oil	274829	1
Caelo	Clove oil	28051313	1
Bombastus	Clove oil	285572	1
Bombastus	Clove oil	288875	2
Bombastus	Clove oil	288127	1
Bombastus	Clove oil	294000	3
Bombastus	Clove oil	295672	2
Bombastus	Clove oil	296428	1
Caelo	Clove oil	2965568	1
Bombastus	Clove oil	296568	3
Bombastus	Clove oil	296979	1
phönix	Clove oil	297263	1
Bombastus	Clove oil	297263	1
Caelo	Clove oil	303337	1
Caelo	Clove oil	30111304	1
Bombastus	Clove oil	300901	1
Bombastus	Clove oil	307074	1
Bombastus	Clove oil	308728	1
Bombastus	Clove oil	308074	1
Bombastus	Clove oil	310830	1
Fagron	Clove oil	3113Q01969	1
Caelo	Clove oil	3512M-01969	1
Caelo	Clove oil	3613Q01969	1
Caelo	Clove oil	402109	1
Taoasis/ Sanacorp	Clove oil	B130612-112414	1
Euro OTC	Clove oil	5210A-01969	1
Taoasis	Clove oil	97204	1
Fagron	Clove oil	826188	1
Wolf	Clove oil	6811	1
Caelo	Clove oil	0-130513-1	1
Caelo	Clove oil	11231705	1
Caelo	Clove oil	12165904	1
Fagron	Clove oil	13L03-B01-292270	1
Caelo	Clove oil	11K17-No7	1
Fagron	Clove oil	11L21-N07	1
Fagron	Clove oil	12C27N04	1
Caelo	Clove oil	13182301-3373319042	1
Bombastus	Clove oil	60010610	1

- 90 918 spectra from 501 *Apo-Ident* customers from a total of 26 879 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Clove oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Clove oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	200	0	93 674
Type C	0	420	14	90 918

The substance/substance group *Clove oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.0000 %)
Type C	100.0000 % (> 99.6586 %)	96.7742 % (> 96.0829 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30858	31707	3.50	41.65
31299	31707	4.54	36.65
31409	31707	5.46	36.07
31526	31526	0.00	27.74
31707	31707	0.00	39.00

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Coal tar spirit SR**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31126-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Coal tar spirit SR; Lithanthracis picis spiritus SR; Spiritus picis lithanthracis SR

Special notes

When selecting the *Coal tar spirit SR* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Coal tar spirit SR*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Coal tar spirit SR	13D05-M14	31126	40	1403343
apomix	Coal tar spirit SR	04C151014	31626	40	20141022
apomix	Coal tar spirit SR	05B170715	33466	40	20150728

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Coal tar spirit SR*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Coal tar spirit SR*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Coal tar spirit SR	2911015-01	30220	40
apomix	Coal tar spirit SR	04C151014	31626 [†]	20
Fagron	Coal tar spirit SR	04B200916	33464	40
Fagron	Coal tar spirit SR	04B200916	33465	40

- 93 734 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 70 spectra from 50 *Apo-Ident* customers from 48 batches from the substance/substance group *Coal tar spirit SR*.
- These include spectra of independent samples from 47 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Gehe	Coal tar spirit SR	2B040413	1
Fagron	Coal tar spirit SR	11E18-M02	1
Apomix /AllianceHC	Coal tar spirit SR	1A140115	1
apomix / Anzag	Coal tar spirit SR	1A140115	1
Apomix	Coal tar spirit SR	1A190216	2
Apomix	Coal tar spirit SR	1B120214	2
Fagron	Coal tar spirit SR	2A070316	4
Fagron	Coal tar spirit SR	2A270215	1
PKH GmbH/Noweda	Coal tar spirit SR	2A070316	1
Caelo	Coal tar spirit SR	2A070316	1
AHD /apomix	Coal tar spirit SR	2B270215	2
Fagron	Coal tar spirit SR	3A130614	1
PKH Halle/Sanacorp	Coal tar spirit SR	3A010415	1
Apomix	Coal tar spirit SR	2C280414	1
pkh halle,phönix	Coal tar spirit SR	4A170615	1
Gehe	Coal tar spirit SR	3B130614	1
Sanacorp	Coal tar spirit SR	4B151014	1
apomix PKH GmbH Halle	Coal tar spirit SR	4A170615	1
AHD /apomix	Coal tar spirit SR	4B200916	1
gehe	Coal tar spirit SR	4B220813	1
Bombastus	Coal tar spirit SR	4C151014	1
PKH Halle / Anzag	Coal tar spirit SR	4C151014	1
Apomix/Noweda	Coal tar spirit SR	4C151014	1
Apomix	Coal tar spirit SR	5B170715	1
PKH Halle	Coal tar spirit SR	4C151014	1
apomix	Coal tar spirit SR	5B170715	1
Apomix / Alliance Healt...	Coal tar spirit SR	5B170715	1
PKH Halle	Coal tar spirit SR	5B170715	1
PKH GmbH	Coal tar spirit SR	5B181113	1
Sanacorp/ apomix	Coal tar spirit SR	5B231116	1
Apomix	Coal tar spirit SR	5B271112	1
AHD /apomix	Coal tar spirit SR	5C231116	1
Apomix	Coal tar spirit SR	5C231116	1
Caelo	Coal tar spirit SR	6A031213	1
Caelo	Coal tar spirit SR	6B031213	1
Apomix Anzag	Coal tar spirit SR	6B031213	1
Fagron	Coal tar spirit SR	6b210915	1
Fagron	Coal tar spirit SR	6B210915	1
Fagron	Coal tar spirit SR	13G16-M01	1
Fagron	Coal tar spirit SR	1404-M02	1
Fagron	Coal tar spirit SR	14I04-M02	1
Fagron	Coal tar spirit SR	14A30-M01	2
Fagron	Coal tar spirit SR	14/04-M02	1
Fagron	Coal tar spirit SR	15C11-M09	1
Fagron	Coal tar spirit SR	170209SKTS02A070316	1
Apomix	Coal tar spirit SR	1A040112	1
Apomix / Sanacorp	Coal tar spirit SR	2C040413	1
apomix	Coal tar spirit SR	3b170412	1
Gehe	Coal tar spirit SR	2B290212	1
Apomix Anzag	Coal tar spirit SR	3B270613	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix	Coal tar spirit SR	3C230508	1
apomix 004	Coal tar spirit SR	4A190712	1
apomix/Anzag	Coal tar spirit SR	4B220813	1
ApoMix	Coal tar spirit SR	4B221009	1
Caelo	Coal tar spirit SR	5C271112	1
Apomix	Coal tar spirit SR	5A101210	1
Apomix/ Noweda	Coal tar spirit SR	6A211212	1
Gehe	Coal tar spirit SR	13D05-M14	1
Gehe	Coal tar spirit SR	6B211212	1
Fagron	Coal tar spirit SR	13D05-M14	1
Caelo	Coal tar spirit SR	15C11-M09	1
Euro OTC	Coal tar spirit SR	1414/04-M02	1
apomix	Coal tar spirit SR	7633599	1

- 91 282 spectra from 501 *Apo-Ident* customers from a total of 26 962 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Coal tar spirit SR* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Coal tar spirit SR* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	103	37	93 734
Type C	0	68	2	91 282

The substance/substance group *Coal tar spirit SR* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	73.5714 % (> 71.4286 %)
Type C	100.0000 % (> 99.6587 %)	97.1429 % (> 92.8571 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31126	31126	0.00	98.21
31626	31126	5.36	98.13
33466	31126	12.38	93.57

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Collodion 4%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30210-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Collodion 4%

Special notes

When selecting the *Collodion 4%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Collodion 4%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Collodion 4%	13325403	31186	40	1402569
Caelo	Collodion 4%	13325411	31434	40	20131016
Caelo	Collodion 4%	16284201	32962	40	20161107

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Collodion 4%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Collodion 4%*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Collodion 4%	09247401	30210	40
Caelo	Collodion 4%	11336004	30533	40
Caelo	Collodion 4%	13325411	31434 [†]	20
Caelo	Collodion 4%	16337201	33514	40
Caelo	Collodion 4%	18070501	34070	40

- 93 694 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 57 spectra from 27 *Apo-Ident* customers from 22 batches from the substance/substance group *Collodion 4%*.
- These include spectra of independent samples from 19 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Collodion 4%	11336004	1
Caelo	Collodion 4%	7031404	1
Caelo	Collodion 4%	13325402	3
Caelo	Collodion 4%	13325403	3
Caelo	Collodion 4%	13325405	3
Caelo	Collodion 4%	13325411	2
Caelo	Collodion 4%	14079101	4
Caelo	Collodion 4%	14079105	3
Euro OTC	Collodion 4%	14079105	1
Caelo	Collodion 4%	15199203	5
Caesar&Loretz GmbH /Anzag	Collodion 4%	14079105	1
Euro OTC	Collodion 4%	15199203	1
Caelo	Collodion 4%	15199206	3
Caelo	Collodion 4%	15429002	1
Caesar & Loretz/Sanacorp	Collodion 4%	15199206	1
Caelo	Collodion 4%	16284201	2
Caelo	Collodion 4%	16337201	3
Caelo	Collodion 4%	27041502	1
Caelo	Collodion 4%	11336001	2
Caelo	Collodion 4%	12035903	2
Caesar&loretz,Phoenix	Collodion 4%	12035916	1
Caelo	Collodion 4%	12220003	1
Defektur	Collodion 4%	12035916	1
Caelo	Collodion 4%	12035913	2
Caelo	Collodion 4%	12035916	5
Caelo	Collodion 4%	90266259	1
Caelo	Collodion 4%	20021305	1
Caelo	Collodion 4%	13325401	1
Gehe	Collodion 4%	11336004	1

- 91 295 spectra from 501 *Apo-Ident* customers from a total of 26 987 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Collodion 4%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Collodion 4%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	167	13	93 694
Type C	0	54	3	91 295

The substance/substance group *Collodion 4%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	92.7778 % (> 91.1111 %)
Type C	100.0000 % (> 99.6587 %)	94.7368 % (> 89.4737 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31186	31186	0.00	194.56
31434	31186	10.99	196.14
32962	31186	27.15	183.72

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Cooling cream / Cooling cream with rose oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30889-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Cooling cream / Cooling cream with rose oil; Cooling cream DAB; Cooling cream DAB 6 (stabilised, contains rose oil); Unguentum leniens; Unguentum leniens DAB 6

Special notes

When selecting the *Cooling cream / Cooling cream with rose oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cooling cream / Cooling cream with rose oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cooling cream DAB	12364802	30889	120	20121217
Bombastus	Cooling cream DAB	275667	30890	80	20121221
Henry Lamotte	Cooling cream DAB	5321500	30921	80	20120627
Caelo	Cooling cream DAB	13051601	31052	40	20130221

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cooling cream DAB	13315101	31224	70	20130925
Caelo	Cooling cream DA...	13317001	31249	40	1405180
Caelo	Cooling cream DA...	14085901	31457	40	20140326
Caelo	Cooling cream DAB	14118002	31471	40	AR-14-FG-012400-01
Caelo	Cooling cream DA...	153029	31922	40	20150911*
Caelo	Cooling cream DAB	153277	31992	40	20150930*
Caelo	Cooling cream DA...	160896	32294	40	20160405*
Caelo	Cooling cream DA...	162223	32506	40	20160727*
Caelo	Cooling cream DAB	172551	33779	30	20170927*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 700 spectra of 13 reference samples from the substance/substance group *Cooling cream / Cooling cream with rose oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 13 different batches.
- 34 844 spectra from a total of 827 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 3072 spectra from 79 reference samples of the substance/substance group *Cooling cream / Cooling cream with rose oil*.
- These include spectra of independent samples from 70 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cooling cream DAB	11107002	30324	42
Caelo	Cooling cream DAB	12001901	30464	30
Infectopharm	Cooling cream DAB 6 (stabili...	12003001	30585	40
Caelo	Cooling cream DAB	12364802	30889 [†]	10
Bombastus	Cooling cream DAB	275667	30890 [†]	10
Caelo	Cooling cream DAB 6 (stabili...	12306001	31037	40
Caelo	Cooling cream DAB 6 (stabili...	153029	31922 [†]	20
Caelo	Cooling cream DAB	153277	31992 [†]	20
Caelo	Cooling cream DAB	153281	31993	60
Caelo	Cooling cream DAB	161293	32227	40
Caelo	Cooling cream DAB	161293	32228	40
Caelo	Cooling cream DAB	161294	32229	40
Caelo	Cooling cream DAB	161295	32230	40
Caelo	Cooling cream DAB	161296	32231	40
Caelo	Cooling cream DAB	161297	32232	40
Caelo	Cooling cream DAB	161341	32233	40
Caelo	Cooling cream DAB	161342	32234	40
Caelo	Cooling cream DAB	161338	32235	40
Caelo	Cooling cream DAB 6 (stabili...	160897	32293	40
Caelo	Cooling cream DAB	161718	32408	40
Caelo	Cooling cream DAB	161719	32409	40
Caelo	Cooling cream DAB	161719	32410	40
Caelo	Cooling cream DAB	161720	32411	40
Caelo	Cooling cream DAB	161720	32412	40
Caelo	Cooling cream DAB	161725	32413	40
Caelo	Cooling cream DAB	161727	32414	40
Caelo	Cooling cream DAB	161728	32415	40
Caelo	Cooling cream DAB	161729	32416	40
Caelo	Cooling cream DAB 6 (stabili...	162224	32507	40
Caelo	Cooling cream DAB	162445	32624	40
Caelo	Cooling cream DAB	162446	32625	40
Caelo	Cooling cream DAB	162447	32626	40
Caelo	Cooling cream DAB	162462	32627	40
Caelo	Cooling cream DAB	162463	32628	40
Caelo	Cooling cream DAB	162464	32629	40
Caelo	Cooling cream DAB	163017	32888	40
Caelo	Cooling cream DAB	163018	32889	40
Caelo	Cooling cream DAB	163019	32890	40
Caelo	Cooling cream DAB	163020	32891	40
Caelo	Cooling cream DAB	163021	32892	40
Caelo	Cooling cream DAB	170004	33112	40
Caelo	Cooling cream DAB	170005	33113	40
Caelo	Cooling cream DAB	163540	33164	40
Caelo	Cooling cream DAB	163544	33165	40
Caelo	Cooling cream DAB	170004	33166	40
Caelo	Cooling cream DAB	170005	33167	40
Caelo	Cooling cream DAB	170003	33168	40
Caelo	Cooling cream DAB	170006	33169	40
Caelo	Cooling cream DAB 6 (stabili...	170413	33180	40
Caelo	Cooling cream DAB	171251	33340	40
Caelo	Cooling cream DAB	171157	33341	40
Caelo	Cooling cream DAB	171250	33342	40
Caelo	Cooling cream DAB	171159	33343	40
Caelo	Cooling cream DAB	171249	33344	40
Caelo	Cooling cream DAB	171160	33345	40
Caelo	Cooling cream DAB	171227	33346	40
Caelo	Cooling cream DAB	171228	33347	40
Caelo	Cooling cream DAB	171156	33348	40
Caelo	Cooling cream DAB	171248	33349	40
Caelo	Cooling cream DAB	170748	33404	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Cooling cream DAB	171567	33563	40
Caelo	Cooling cream DAB	171570	33564	40
Caelo	Cooling cream DAB	171571	33565	40
Caelo	Cooling cream DAB	171572	33566	40
Caelo	Cooling cream DAB	171573	33567	40
Caelo	Cooling cream DAB	171897	33568	40
Caelo	Cooling cream DAB	172536	33771	40
Caelo	Cooling cream DAB	172539	33773	40
Caelo	Cooling cream DAB	172540	33774	40
Caelo	Cooling cream DAB	172541	33775	40
Caelo	Cooling cream DAB	172544	33776	40
Caelo	Cooling cream DAB	172545	33777	40
Caelo	Cooling cream DAB	172546	33778	40
Caelo	Cooling cream DAB	172548	33780	40
Caelo	Cooling cream DAB	172549	33781	40
Caelo	Cooling cream DAB	172550	33782	40
Caelo	Cooling cream DAB 6 (stabili...	180218	34142	40
Caelo	Cooling cream DAB	180424	34150	40
Caelo	Cooling cream DAB 6 (stabili...	181897	34265	40

- 90 802 spectra from a total of 2093 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 3425 spectra from 383 *Apo-Ident* customers from 988 batches from the substance/substance group *Cooling cream / Cooling cream with rose oil*.
- These include spectra of independent samples from 976 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	12011402	1
Caelo	Cooling cream DAB	12090401	3
Caelo	Cooling cream DAB	12292401	4
Caelo	Cooling cream DAB	12364601	14
Caelo	Cooling cream DAB	12312102	5
Caelo	Cooling cream DAB	13051601	9
Caelo	Cooling cream DAB	13059901	12
Phönix	Cooling cream DAB	13121001	3
Noweda	Cooling cream DAB	13121002	1
Phönix	Cooling cream DAB	13121701	1
Caelo	Cooling cream DAB	13175101	7
Caelo	Cooling cream DAB	13175201	7

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	13191501	4
Caelo	Cooling cream DAB	13191701	7
Phönix	Cooling cream DAB	13191801	1
Caelo	Cooling cream DAB	13191901	4
Caelo	Cooling cream DAB	13192101	9
Caelo	Cooling cream DAB	13191801	8
Berkel AHK	Cooling cream DAB	13226301	3
Caelo	Cooling cream DAB	13226501	9
Caelo	Cooling cream DAB	13226801	23
gehe 7.12.13	Cooling cream DAB	13226701	1
Phönix	Cooling cream DAB	13227401	1
Caelo	Cooling cream DAB	13315701	3
Caelo	Cooling cream DAB	13315101	28
Caelo	Cooling cream DAB	13402201	16
Caelo	Cooling cream DAB	14000201	7
Caelo	Cooling cream DAB	14000301	12
Phönix	Cooling cream DAB	14055401	2
Caelo	Cooling cream DAB	14000401	2
Caelo	Cooling cream DAB	14055301	7
Caelo	Cooling cream DAB	14086201	5
Caelo	Cooling cream DAB	14129401	14
Noweda	Cooling cream DAB	14085901	1
Großhandel	Cooling cream DAB	14055501	1
Caelo	Cooling cream DAB	14086301	13
Caelo	Cooling cream DAB	14129801	20
Caelo	Cooling cream DAB	14145701	5
Caelo	Cooling cream DAB	14145901	10
Caelo	Cooling cream DAB	14184501	17
Caelo	Cooling cream DAB	14184601	25
Caelo	Cooling cream DAB	14185101	5
Caelo	Cooling cream DAB	14268101	7
Caelo	Cooling cream DAB	142681012	1
Caelo	Cooling cream DAB	14278201	11
Caelo	Cooling cream DAB	14271101	4
Caelo	Cooling cream DAB	14278301	17
Caelo	Cooling cream DAB	14278502	5
Caelo	Cooling cream DAB	14328702	8
Caelo	Cooling cream DAB	14328801	18
Caelo	Cooling cream DAB	15052801	19
Caelo	Cooling cream DAB	15153701	8
Caelo	Cooling cream DAB	15228601	35
Caelo	Cooling cream DAB	15152801	3
Caelo	Cooling cream DAB	16000201	6
Caelo	Cooling cream DAB	15435101	10
Caelo	Cooling cream DAB	16000801	4
Caelo	Cooling cream DAB	1915M-02541	1
Caelo	Cooling cream DAB	16129501	8
Bombastus	Cooling cream DAB	275666	2
Caelo	Cooling cream DAB	279753	1
Ginkgo-Apotheke	Cooling cream DAB	278439	1
Caelo	Cooling cream DAB	16172601	19
Henry Lamotte/GEHE	Cooling cream DAB	5407300	1
lamotte/sanacorb	Cooling cream DAB	5496200	1
Henry Lamotte Oils/GEHE	Cooling cream DAB	5496200	1
Lamotte/AHD	Cooling cream DAB	5600300	2
Sanacorp	Cooling cream DAB		1
sanacorp/Lamotte	Cooling cream DAB	5735600	1
Sanacorp	Cooling cream DAB	76/AKD3	1
Caelo	Cooling cream DAB	73/AKD3	1
Sanacorp	Cooling cream DAB	73AKD4	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	14184801	5
Bombastus	Cooling cream DAB 6 (stabili...	2/2018	1
Caelo	Cooling cream DAB	123456	1
Bombastus	Cooling cream DAB	104948	1
Caelo	Cooling cream DAB	123test	2
Caelo	Cooling cream DAB	1256856	1
Caelo	Cooling cream DAB 6 (stabili...	15328001	8
Caelo	Cooling cream DAB	15434402	11
Caelo	Cooling cream DAB	15434502	8
Caelo	Cooling cream DAB	15434501	14
Caelo	Cooling cream DAB 6 (stabili...	15435101	2
Caelo	Cooling cream DAB	15435102	16
Caelo	Cooling cream DAB 6 (stabili...	16000201	3
Caelo	Cooling cream DAB	16000501	5
Caelo	Cooling cream DAB	16000601	11
Phoenix	Cooling cream DAB 6 (stabili...	16000501	2
Ichthyol	Cooling cream DAB	16000801	1
Phoenix	Cooling cream DAB	16085801	1
Caelo	Cooling cream DAB 6 (stabili...	16085801	1
Caelo	Cooling cream DAB	1605901	1
Caelo	Cooling cream DAB	16085801	6
Caelo	Cooling cream DAB	16085901	6
Caelo	Cooling cream DAB 6 (stabili...	16085901	3
Caelo	Cooling cream DAB	16085902	4
Caelo	Cooling cream DAB	16086001	5
sanacorp/ichthyol	Cooling cream DAB	16086001	1
Phoenix	Cooling cream DAB	16086001	1
Caelo	Cooling cream DAB	16086101	7
Caelo	Cooling cream DAB 6 (stabili...	16086001	1
Caelo	Cooling cream DAB 6 (stabili...	16086101	1
Bombastus	Cooling cream DAB	16095503	1
Caelo	Cooling cream DAB	16129401	7
Caelo	Cooling cream DAB 6 (stabili...	16129501	3
Euro OTC	Cooling cream DAB	16129501	1
Caelo	Cooling cream DAB	16129502	1
Caelo	Cooling cream DAB	16129701	3
Caelo	Cooling cream DAB 6 (stabili...	16129701	2
Caelo	Cooling cream DAB	16133901	3
Caelo	Cooling cream DAB	16134001	3
Phoenix	Cooling cream DAB	16134001	1
Caelo	Cooling cream DAB 6 (stabili...	16134001	1
Caelo	Cooling cream DAB	16134101	3
Caelo	Cooling cream DAB 6 (stabili...	16134101	3
Caelo	Cooling cream DAB	16134201	17
Euro OTC	Cooling cream DAB	16134201	1
Caelo	Cooling cream DAB 6 (stabili...	16134201	2
Caelo	Cooling cream DAB	16171901	4
Caesar & Loretz/Noweda	Cooling cream DAB 6 (stabili...	16134201	1
Caelo	Cooling cream DAB 6 (stabili...	16171901	1
Phoenix	Cooling cream DAB	16172001	1
Caelo	Cooling cream DAB	16172001	6
Caelo	Cooling cream DAB 6 (stabili...	16172001	1
Celo/ Kehr	Cooling cream DAB	16172001	1
Caelo	Cooling cream DAB 6 (stabili...	16172501	1
Ichthyol/Phönix	Cooling cream DAB 6 (stabili...	16172601	2
Ichthyol-Gesellschaft	Cooling cream DAB	16172601	1
Caelo	Cooling cream DAB 6 (stabili...	16172601	9
Caelo	Cooling cream DAB	16172701	4
Caelo	Cooling cream DAB 6 (stabili...	16172701	2
Phoenix	Cooling cream DAB	16172701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	16172801	5
Caelo	Cooling cream DAB	16172901	8
Caelo	Cooling cream DAB 6 (stabili...	16172901	4
Caelo	Cooling cream DAB	16172902	5
Caelo	Cooling cream DAB 6 (stabili...	16172902	1
Caelo	Cooling cream DAB 6 (stabili...	16244601	1
Gehe	Cooling cream DAB	16244601	1
Caelo	Cooling cream DAB	16244601	6
Caelo	Cooling cream DAB 6 (stabili...	16244701	9
Caelo	Cooling cream DAB	16244701	31
Caelo	Cooling cream DAB	16246201	4
Caelo	Cooling cream DAB	16246301	5
Caelo	Cooling cream DAB 6 (stabili...	16246301	2
Caelo	Cooling cream DAB	16301702	3
Caelo	Cooling cream DAB	16301801	17
Caelo	Cooling cream DAB 6 (stabili...	16301702	1
Sanacorp 08.03.2017 EK:...	Cooling cream DAB	16301801	1
Fagron	Cooling cream DAB	16301801	1
Caelo	Cooling cream DAB	16301901	3
Gehe	Cooling cream DAB 6 (stabili...	16301901	1
Caelo	Cooling cream DAB	16302001	4
Gehe	Cooling cream DAB 6 (stabili...	16302001	1
Caelo	Cooling cream DAB 6 (stabili...	16302001	1
Caelo	Cooling cream DAB	16302101	2
Caelo	Cooling cream DAB	16322901	2
Caelo	Cooling cream DAB	163230001	1
Caelo	Cooling cream DAB	16323001	5
Noweda	Cooling cream DAB	16323101	1
Caelo	Cooling cream DAB	16323101	3
Caesar&Loretz	Cooling cream DAB	16354401	2
Caelo/Phönix	Cooling cream DAB	16354401	1
Caelo	Cooling cream DAB	16354401	9
AHD	Cooling cream DAB	16354401	2
Caelo	Cooling cream DAB 6 (stabili...	16354401	3
Caelo	Cooling cream DAB	16354402	4
Caelo	Cooling cream DAB	17000301	4
Caelo	Cooling cream DAB	17000601	24
Caelo	Cooling cream DAB	17000401	2
Caelo	Cooling cream DAB 6 (stabili...	17000401	1
Bombastus	Cooling cream DAB	17000979	1
Bombastus	Cooling cream DAB	17000982	2
Bombastus	Cooling cream DAB	17000976	1
Finteler Apotheke	Cooling cream DAB 6 (stabili...	17000601	1
Bombastus	Cooling cream DAB	17000983	1
Bombastus	Cooling cream DAB	17001962	1
Bombastus	Cooling cream DAB	17001964	1
Caelo	Cooling cream DAB	170207UL/16323101	1
Bombastus	Cooling cream DAB	17001963	1
Caelo	Cooling cream DAB	17033001	15
hochhaus-Apotheke	Cooling cream DAB 6 (stabili...	17033001	1
Caelo	Cooling cream DAB 6 (stabili...	17033001	4
Caelo	Cooling cream DAB 6 (stabili...	17035301	2
Caelo	Cooling cream DAB	17035301	7
Caelo	Cooling cream DAB	17035401	2
Caelo	Cooling cream DAB 6 (stabili...	17035401	1
Caelo	Cooling cream DAB	17038501	7
Gehe	Cooling cream DAB 6 (stabili...	17035501	1
Phoenix	Cooling cream DAB	17038501	1
Caelo	Cooling cream DAB	170728UL/17081801	1
Caelo	Cooling cream DAB 6 (stabili...	17074801	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB 6 (stabili...	17075001	15
Caelo	Cooling cream DAB	17075101	5
Phönix	Cooling cream DAB 6 (stabili...	17075101	1
Caelo	Cooling cream DAB 6 (stabili...	17075101	3
Noweda	Cooling cream DAB	17075101	1
Caelo	Cooling cream DAB	17075201	10
phönix	Cooling cream DAB 6 (stabili...	17075201	1
Caelo	Cooling cream DAB 6 (stabili...	17075201	3
Caelo	Cooling cream DAB	17075301	6
Caelo	Cooling cream DAB 6 (stabili...	17075301	3
Caelo	Cooling cream DAB	17075401	4
Caelo	Cooling cream DAB 6 (stabili...	17081801	3
Caelo	Cooling cream DAB	17081801	4
Caelo	Cooling cream DAB	170926UL/17122701	1
Caelo	Cooling cream DAB	17115601	6
Caelo	Cooling cream DAB 6 (stabili...	17115601	3
Caelo	Cooling cream DAB 6 (stabili...	17115701	13
Caelo	Cooling cream DAB 6 (stabili...	17115901	1
Caelo	Cooling cream DAB	17115901	1
Caelo	Cooling cream DAB	17122701	2
Caelo	Cooling cream DAB	17116001	1
Noweda	Cooling cream DAB	17122801	2
Caesar und Loretz GmbH	Cooling cream DAB	17122801	1
Caelo	Cooling cream DAB 6 (stabili...	17124801	4
Caelo	Cooling cream DAB	171229UL/17122801	1
Caelo	Cooling cream DAB 6 (stabili...	1715002	1
Caelo	Cooling cream DAB	17156901	1
Caelo	Cooling cream DAB	17125002	3
Sanacorp	Cooling cream DAB	17157001	1
Caelo	Cooling cream DAB	17157001	2
Caesar und Loretz GmbH	Cooling cream DAB 6 (stabili...	17157201	1
Caelo	Cooling cream DAB	17157101	2
Sanacorp W.E: 18.01.18 ...	Cooling cream DAB 6 (stabili...	17157001	2
Caelo	Cooling cream DAB	17157301	2
Caelo	Cooling cream DAB 6 (stabili...	17157301	1
Alliance	Cooling cream DAB	17189701	1
Caelo	Cooling cream DAB	17189701	3
Caelo	Cooling cream DAB	17189802	1
AHD	Cooling cream DAB	17189701	1
chaelo/gehe	Cooling cream DAB	17254101	1
GEHE	Cooling cream DAB	17254101	1
Caelo	Cooling cream DAB	17254201	10
Caelo	Cooling cream DAB	17254401	1
Caelo	Cooling cream DAB	17254601	1
Caelo	Cooling cream DAB	17298001	6
Caelo	Cooling cream DAB	17254801	1
Caelo	Cooling cream DAB	17255001	1
Caelo	Cooling cream DAB 6 (stabili...	17254501	1
Caelo	Cooling cream DAB	17255101	2
Caelo	Cooling cream DAB	17298101	2
Caelo	Cooling cream DAB	17298401	1
Caelo	Cooling cream DAB		5
Caelo	Cooling cream DAB	17298501	2
Caelo	Cooling cream DAB 6 (stabili...	17298701	1
Caelo	Cooling cream DAB	17330801	1
Caelo	Cooling cream DAB	17330802	4
Caelo	Cooling cream DAB 6 (stabili...	17331001	1
Bombastus	Cooling cream DAB	18001078	1
Bombastus	Cooling cream DAB	18000333	1
Bombastus	Cooling cream DAB	18002087	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB 6 (stabili...	18021803	1
Caelo	Cooling cream DAB	18042701	1
Bombastus	Cooling cream DAB	299351	1
Bombastus	Cooling cream DAB	298117	1
Bombastus	Cooling cream DAB	299362	1
Bombastus	Cooling cream DAB	299967	1
Bombastus	Cooling cream DAB	299968	1
Euro OTC	Cooling cream DAB	299976	1
Bombastus	Cooling cream DAB	299980	2
Bombastus	Cooling cream DAB	299981	1
Bombastus	Cooling cream DAB	300720	2
Bombastus	Cooling cream DAB	300715	1
Bombastus	Cooling cream DAB	300717	1
Bombastus	Cooling cream DAB	300727	1
Bambustus/Noweda	Cooling cream DAB	300727	2
Bombastus	Cooling cream DAB	300728	2
Bombastus	Cooling cream DAB	300730	2
Bombastus	Cooling cream DAB	301670	1
Bombastus	Cooling cream DAB	301672	1
Bombastus	Cooling cream DAB	301674	1
Bombastus	Cooling cream DAB	301675	2
Bombastus	Cooling cream DAB	301699	2
Bombastus	Cooling cream DAB	301700	3
Bombastus	Cooling cream DAB	301701	3
Ebert & Jacobi 02.08.16	Cooling cream DAB	301701	1
Fagron	Cooling cream DAB	302370	2
Caelo	Cooling cream DAB	301701	1
Bombastus	Cooling cream DAB	302371	1
Bombastus	Cooling cream DAB	302372	1
Bombastus	Cooling cream DAB	302374	2
Bombastus	Cooling cream DAB	302380	1
Bombastus	Cooling cream DAB	302375	1
Bombastus	Cooling cream DAB	302385	4
Bombastus	Cooling cream DAB	303153	1
GEHE	Cooling cream DAB	303163	2
Bombastus	Cooling cream DAB 6 (stabili...	303620	1
Bombastus	Cooling cream DAB	303620	1
Bombastus	Cooling cream DAB	303621	2
Bombastus	Cooling cream DAB 6 (stabili...	303621	1
Bombastus	Cooling cream DAB	304582	2
Bombastus	Cooling cream DAB	303625	1
Bombastus	Cooling cream DAB 6 (stabili...	304585	1
Caelo	Cooling cream DAB	304582	1
Bombastus	Cooling cream DAB	304585	1
Bombastus/Phönix	Cooling cream DAB	304582	1
Caelo	Cooling cream DAB	304590	1
Bombastus	Cooling cream DAB	304591	1
Caelo	Cooling cream DAB	304591	1
Bombastus	Cooling cream DAB	304593	3
GEHE	Cooling cream DAB	304594	1
Bombastus	Cooling cream DAB	305017	2
Bombastus	Cooling cream DAB	305021	1
Fagron	Cooling cream DAB	305129	2
Bombastus	Cooling cream DAB	305022	1
Bombastus	Cooling cream DAB	305131	3
Bombastus	Cooling cream DAB	305952	1
Bombastus	Cooling cream DAB	305953	1
GEHE	Cooling cream DAB	305953	1
Bombastus	Cooling cream DAB	305961	2
Bombastus	Cooling cream DAB	305954	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Cooling cream DAB	306529	3
Bombastus	Cooling cream DAB	306530	1
Bombastus	Cooling cream DAB	306532	1
Euro OTC	Cooling cream DAB	306931	1
Bombastus	Cooling cream DAB	306932	3
Phoenix	Cooling cream DAB	306932	1
Caelo	Cooling cream DAB	306941	1
Bombastus	Cooling cream DAB	306940	1
Bombastus	Cooling cream DAB	306944	1
Bombastuts/Gehe	Cooling cream DAB	307350	1
Caelo	Cooling cream DAB	307356	1
Bombastus	Cooling cream DAB	307916	2
Bombastus	Cooling cream DAB	307358	1
Bombastuts/Gehe	Cooling cream DAB	307919	1
Bombastus	Cooling cream DAB	307929	3
Bombastus	Cooling cream DAB	307931	1
Bombastus	Cooling cream DAB	307932	1
Bombastus	Cooling cream DAB	308489	1
GEHE	Cooling cream DAB	307933	1
Bombastus	Cooling cream DAB	308491	1
Bombastus	Cooling cream DAB	308710	1
Bombastus	Cooling cream DAB 6 (stabili...	308711	1
Bombastuts/Gehe	Cooling cream DAB	308712	1
Bombastus	Cooling cream DAB	308715	1
GEHE	Cooling cream DAB	308725	1
GEHE	Cooling cream DAB	309370	1
Bombastus	Cooling cream DAB	309371	1
Bombastus	Cooling cream DAB	310000	2
Bombastus	Cooling cream DAB	310002	1
Bombastus	Cooling cream DAB	310974	3
Bombastus	Cooling cream DAB	310977	1
Bombastus	Cooling cream DAB	310975	1
Bombastus	Cooling cream DAB	310984	1
Lamotte/Noweda	Cooling cream DAB	6002062003	1
Sanacorp/Henry Lamotte	Cooling cream DAB	6003682010	3
Bombastus	Cooling cream DAB	310987	1
Henry Lamtte Oils GmbH	Cooling cream DAB	6002062002	1
Bombastus	Cooling cream DAB	60082801	2
Lamotte/Noweda	Cooling cream DAB 6 (stabili...	8002112004	1
Lamotte/Noweda	Cooling cream DAB	8002320008	1
Henry Lamotte/Anzag	Cooling cream DAB	8002320007	1
Caelo	Cooling cream DAB	8002320007	2
Ebert/Lamotte	Cooling cream DAB	8002696003	1
Gehe	Cooling cream DAB	8002320008	1
Henry Lamotte	Cooling cream DAB	8003158004	1
Henry Lamotte Oils / Gehe	Cooling cream DAB	8002875004	1
Henry Lamotte/Noweda	Cooling cream DAB	8003478005	1
Caelo	Cooling cream DAB	8002875004	1
Caelo	Cooling cream DAB	1097065	1
Caelo	Cooling cream DAB 6 (stabili...	1097214	1
Caelo	Cooling cream DAB 6 (stabili...	15228501	4
WE:24.02.2016 Phönix ...	Cooling cream DAB	15228501	2
Caelo	Cooling cream DAB 6 (stabili...	15237101	4
Caelo	Cooling cream DAB	15228501	15
Caelo	Cooling cream DAB	15328001	31
Caesar & Loretz/Noweda	Cooling cream DAB 6 (stabili...	15328001	2
Caelo	Cooling cream DAB 6 (stabili...	15417301	2
Caelo	Cooling cream DAB	15413901	3
Caelo	Cooling cream DAB	15417401	2
Caelo	Cooling cream DAB	15417501	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alliance H.	Cooling cream DAB	15434401	1
Caelo	Cooling cream DAB	15434401	15
Caelo	Cooling cream DAB 6 (stabili...	15434501	7
Caelo	Cooling cream DAB 6 (stabili...	15434502	6
Phönix	Cooling cream DAB 6 (stabili...	15434502	1
Caelo	Cooling cream DAB 6 (stabili...	15434503	1
Caelo	Cooling cream DAB	15434503	4
Caelo	Cooling cream DAB	1543501	1
Caelo	Cooling cream DAB	1543502	1
Euro OTC	Cooling cream DAB	15435101	1
Caelo	Cooling cream DAB 6 (stabili...	15435102	1
	Cooling cream DAB	15L03-T06-028365	1
Caelo	Cooling cream DAB 6 (stabili...	16000501	1
Noweda	Cooling cream DAB	16000801	1
Caelo	Cooling cream DAB 6 (stabili...	16000601	1
Caelo	Cooling cream DAB 6 (stabili...	16089702	1
Caelo	Cooling cream DAB	16129601	3
Henry Lamotte/ Henry La...	Cooling cream DAB 6 (stabili...	16171901	1
Caesar & Loretz/Noweda	Cooling cream DAB 6 (stabili...	16172601	1
Caelo	Cooling cream DAB 6 (stabili...	16172801	1
Caesar & Loretz GmbH	Cooling cream DAB 6 (stabili...	16222401	1
Caelo	Cooling cream DAB 6 (stabili...	16222401	1
Caelo	Cooling cream DAB 6 (stabili...	16222402	1
Caelo	Cooling cream DAB	162446901	1
Caelo	Cooling cream DAB 6 (stabili...	16301801	2
Caelo	Cooling cream DAB	163020001	1
Bombastus	Cooling cream DAB	17000227	1
Caelo	Cooling cream DAB	17000501	1
Bombastus	Cooling cream DAB	17000973	1
Bombastus	Cooling cream DAB	17001966	1
Gehe	Cooling cream DAB	17002971	1
Caelo	Cooling cream DAB 6 (stabili...	17041301	4
Caelo	Cooling cream DAB 6 (stabili...	170522UL/16323101	1
Caelo	Cooling cream DAB 6 (stabili...	17116001	1
Caelo	Cooling cream DAB 6 (stabili...	17117402	4
Caelo	Cooling cream DAB	17122801	1
Caelo	Cooling cream DAB 6 (stabili...	17125002	1
Caelo	Cooling cream DAB	17156801	1
Caelo	Cooling cream DAB 6 (stabili...	17156901	1
Caelo	Cooling cream DAB 6 (stabili...	17156801	1
Caelo	Cooling cream DAB	17254501	1
Caelo	Cooling cream DAB 6 (stabili...	18106401	2
Caelo	Cooling cream DAB 6 (stabili...	18002802	1
Bombastus	Cooling cream DAB	298119	2
Bombastus	Cooling cream DAB	298113	3
Caelo	Cooling cream DAB	299352	1
Bombastus	Cooling cream DAB	299966	1
Bombastus	Cooling cream DAB	299978	1
Bombastuts/Noweda	Cooling cream DAB	300726	1
Bombastus	Cooling cream DAB	303158	1
Bombastus	Cooling cream DAB	305129	1
Bombastus	Cooling cream DAB	302382	1
Bombastus	Cooling cream DAB	305949	1
Bombastus	Cooling cream DAB	309999	1
Caelo	Cooling cream DAB	5215A-02541	1
Euro OTC	Cooling cream DAB	308491	1
Caelo	Cooling cream DAB	6002062002	1
Henry Lamotte Oils GmbH...	Cooling cream DAB	6002062002	1
lamotte/anzag	Cooling cream DAB	8002112004	1
Bombastus	Cooling cream DAB	310976	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils GmbH...	Cooling cream DAB	8002320007	1
Bombastus	Cooling cream DAB	265549	1
Bombastus	Cooling cream DAB	266717	1
Bombastus	Cooling cream DAB	265548	1
lamotte/anzag	Cooling cream DAB	8002320008	1
Bombastus	Cooling cream DAB	265554	1
Bombastus	Cooling cream DAB	264601	1
Gehe	Cooling cream DAB	266839	2
Bombastus	Cooling cream DAB	266941	1
Bombastus	Cooling cream DAB	268073	3
Bombastus	Cooling cream DAB	268194	2
Phönix	Cooling cream DAB	268978	1
Bombastus	Cooling cream DAB	268981	1
Bombastus	Cooling cream DAB	269088	1
Gehe	Cooling cream DAB	269093	2
Gehe	Cooling cream DAB	269098	3
Gehe	Cooling cream DAB	271535	1
Bombastus	Cooling cream DAB	270381	1
Bombastus	Cooling cream DAB	270309	1
Bombastus	Cooling cream DAB	271632	1
Caelo	Cooling cream DAB	1071407	1
Caelo	Cooling cream DAB	2091404	1
Caelo	Cooling cream DAB	1097071	1
Caelo	Cooling cream DAB	1071408	1
Bombastus	Cooling cream DAB	3022015B	1
Bombastus	Cooling cream DAB	3092014A	1
Bombastus	Cooling cream DAB	4032015B	1
Caelo	Cooling cream DAB	4051503	2
Caelo	Cooling cream DAB	4061203	1
Caelo	Cooling cream DAB	5031509	1
Caelo	Cooling cream DAB	4061405	1
Caelo	Cooling cream DAB	7071203	1
Caelo	Cooling cream DAB	8121403	1
Bombastus	Cooling cream DAB	6861938	1
Caelo	Cooling cream DAB	11268501	1
Caesar & Loretz GmbH 40...	Cooling cream DAB	11299001	1
Caelo	Cooling cream DAB	11321402	1
Caelo	Cooling cream DAB	11322901	1
Caelo	Cooling cream DAB	11368101	1
Caelo	Cooling cream DAB	11323101	2
Caelo	Cooling cream DAB	12002101	2
Caelo	Cooling cream DAB	12001901	3
Caelo	Cooling cream DAB	12002401	1
Caelo	Cooling cream DAB	11368201	1
Caelo	Cooling cream DAB	12002501	3
Fiebig ek: 9.07EUR	Cooling cream DAB	1204001	1
Caelo	Cooling cream DAB	12045901	1
Caesar & Loretz GmbH, 4...	Cooling cream DAB	12046001	1
Noweda	Cooling cream DAB	12045901	1
Noweda	Cooling cream DAB	12046101	2
Noweda	Cooling cream DAB	12046201	1
Caesar & Loretz GmbH, 4...	Cooling cream DAB	12046201	2
Kehr	Cooling cream DAB	12046101	1
Caelo	Cooling cream DAB	12046401	1
Caelo	Cooling cream DAB	12051402	1
Caelo	Cooling cream DAB	12051401	1
Caelo	Cooling cream DAB	12051701	2
Caelo	Cooling cream DAB	120601202	1
DAB 2008	Cooling cream DAB	12089601	1
Caelo	Cooling cream DAB	12089601	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	12089701	1
Caelo	Cooling cream DAB	12089901	1
Kehr	Cooling cream DAB	12089801	1
Anzag 31.10.12 ek:9,12EUR	Cooling cream DAB	12187301	1
Anzag 27.12.12 ek:9,09	Cooling cream DAB	12312301	1
Anzag 07.05.13 ek: 9,38EUR	Cooling cream DAB	13059701	1
Caelo	Cooling cream DAB	13120901	18
Caelo	Cooling cream DAB	1315601	1
Anzag 05.02.13 ek: 9,38EUR	Cooling cream DAB	12312801	1
Caelo	Cooling cream DAB	13175301	24
Caelo	Cooling cream DAB	13214010	1
Caelo	Cooling cream DAB	13226601	29
Sanacorp	Cooling cream DAB	13227001	1
Caelo	Cooling cream DAB 6 (stabili...	1328702	1
Caelo	Cooling cream DAB	13315201	10
Caelo	Cooling cream DAB	13315801	6
Caelo	Cooling cream DAB	13315601	6
Caelo	Cooling cream DAB	13316001	1
Caelo	Cooling cream DAB	13401801	31
Caelo	Cooling cream DAB	13401901	14
Caelo	Cooling cream DAB	13402101	7
Caelo	Cooling cream DAB	14000101	2
Caelo	Cooling cream DAB	14000501	5
Großhandel	Cooling cream DAB	14055301	1
Caelo	Cooling cream DAB	14028501	1
Bombastus	Cooling cream DAB	14052014B	1
Caelo	Cooling cream DAB	14055401	1
Großhandel	Cooling cream DAB	14055401	1
Caelo	Cooling cream DAB	14055501	2
Caelo	Cooling cream DAB	14055601	10
Caelo	Cooling cream DAB	14055702	11
GEHE 14.04.2014	Cooling cream DAB	14055703	2
Caelo	Cooling cream DAB	14071204	1
Bombastus	Cooling cream DAB	14062101	1
Caelo	Cooling cream DAB	14085901	14
Caesar & Loretz GmbH	Cooling cream DAB	14085901	3
Gehe	Cooling cream DAB	14085901	1
Caelo	Cooling cream DAB	14085902	4
Caelo	Cooling cream DAB	14086001	3
Phoenix 20.06.2014	Cooling cream DAB	14086001	1
Caelo	Cooling cream DAB	14086101	3
Caelo	Cooling cream DAB	14086201	1
Caelo	Cooling cream DAB	14086401	3
Caelo	Cooling cream DAB	1408631	1
Caelo	Cooling cream DAB	14086501	1
Caelo	Cooling cream DAB	14086601	3
Caelo	Cooling cream DAB	14086801	4
Caelo	Cooling cream DAB	14101401	1
Caelo	Cooling cream DAB	1412801	1
Caelo	Cooling cream DAB	14097214	1
Gehe	Cooling cream DAB	14129401	2
Gehe	Cooling cream DAB	14129402	1
Caelo	Cooling cream DAB	14129402	2
Caelo	Cooling cream DAB	14129501	4
Noweda	Cooling cream DAB	14129801	2
Fagron	Cooling cream DAB	14132001	1
Caelo	Cooling cream DAB	14132001	1
Caelo	Cooling cream DAB	14132201	4
Caelo	Cooling cream DAB	14132301	8
Caelo	Cooling cream DAB	14132401	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	14143101	1
Caelo	Cooling cream DAB	14145501	25
Sanacorp	Cooling cream DAB	14145501	1
Caelo	Cooling cream DAB	14145602	1
Caelo	Cooling cream DAB	14145801	4
Caelo	Cooling cream DAB	1418401	2
Caelo	Cooling cream DAB	14184401	8
Caelo	Cooling cream DAB	14184301	1
Caelo	Cooling cream DAB 6 (stabili...	14184401	1
Caelo	Cooling cream DAB	14184402	3
Caesar & Loretz GmbH	Cooling cream DAB	14184501	1
Caelo	Cooling cream DAB	14184701	5
Caelo	Cooling cream DAB	14184901	5
Caelo	Cooling cream DAB	14185001	8
Caelo	Cooling cream DAB 6 (stabili...	14184801	1
Caelo	Cooling cream DAB	14200401	1
Caelo	Cooling cream DAB	14258301	1
Caelo	Cooling cream DAB	14186601	1
Caelo	Cooling cream DAB 6 (stabili...	14268101	2
Sanacorp	Cooling cream DAB	14268101	1
Berkel AHK	Cooling cream DAB	14268101	1
Caelo	Cooling cream DAB	14268102	1
Caelo	Cooling cream DAB	14268201	3
Caelo	Cooling cream DAB	14268301	16
Sana 29.12.14 EK 5,59	Cooling cream DAB	14268301	1
DAB 2010	Cooling cream DAB	14269301	1
Allianz Healthcare 05.1...	Cooling cream DAB	14271101	1
Caelo	Cooling cream DAB	14271201	2
Caelo	Cooling cream DAB 6 (stabili...	14271201	1
Caelo	Cooling cream DAB	14271301	2
Caelo	Cooling cream DAB	14271401	7
Caelo	Cooling cream DAB 6 (stabili...	14271401	2
Caelo	Cooling cream DAB	14271501	3
Caelo	Cooling cream DAB	142728301	1
Caelo	Cooling cream DAB	1427801	1
Caelo	Cooling cream DAB	14278401	3
Caelo	Cooling cream DAB 6 (stabili...	14278401	1
Caelo	Cooling cream DAB	14278601	4
Caelo	Cooling cream DAB	14288601	6
Caesar & Loretz GmbH	Cooling cream DAB	14278701	1
Caelo	Cooling cream DAB 6 (stabili...	14288601	1
Caelo	Cooling cream DAB	14288701	4
Caelo	Cooling cream DAB	14288701	1
Caelo/AHD	Cooling cream DAB	14328702	1
Caelo	Cooling cream DAB	14328701	2
Henry Lamotte	Cooling cream DAB	14328702	1
Caelo	Cooling cream DAB 6 (stabili...	14328702	3
Caelo	Cooling cream DAB 6 (stabili...	14328801	4
Sanacorp	Cooling cream DAB	14328801	2
Phönix	Cooling cream DAB	14328801	1
Caelo	Cooling cream DAB	14328901	3
Caelo	Cooling cream DAB	14329001	5
Caelo	Cooling cream DAB 6 (stabili...	14329001	1
Caelo	Cooling cream DAB	14329201	4
Caelo	Cooling cream DAB 6 (stabili...	14329201	1
Caelo	Cooling cream DAB	14329301	7
Caelo	Cooling cream DAB	14329401	8
Caelo	Cooling cream DAB 6 (stabili...	14329401	1
phönix	Cooling cream DAB 6 (stabili...	14329501	1
Gehe	Cooling cream DAB	14329501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	14329501	3
Caelo	Cooling cream DAB	14329601	8
Caelo	Cooling cream DAB 6 (stabili...	14329601	1
Caelo	Cooling cream DAB	1468201	1
Caelo	Cooling cream DAB	1484701	1
Phönix	Cooling cream DAB	14850901	1
Kremer/Phönix	Cooling cream DAB	15/158/K	1
Fagron	Cooling cream DAB	14A20	1
Cealo/ Noweda	Cooling cream DAB	14c06-t02010890	1
Caelo	Cooling cream DAB	15017601	1
Caelo	Cooling cream DAB	1502801	1
Caelo	Cooling cream DAB 6 (stabili...	15014701	1
Caelo	Cooling cream DAB	15049201	20
Caelo	Cooling cream DAB	15048701	1
Caelo	Cooling cream DAB	15049202	2
Caelo	Cooling cream DAB 6 (stabili...	15049202	5
Caelo	Cooling cream DAB 6 (stabili...	15049201	1
Caelo	Cooling cream DAB 6 (stabili...	15052801	2
Caelo	Cooling cream DAB	15053001	3
Caelo	Cooling cream DAB	15058801	4
Caelo	Cooling cream DAB	15058901	3
Caelo	Cooling cream DAB	15059001	3
Caelo	Cooling cream DAB 6 (stabili...	15059001	2
Caelo	Cooling cream DAB	15059201	26
Caelo	Cooling cream DAB 6 (stabili...	15059201	2
Caelo	Cooling cream DAB	15066501	2
Caelo	Cooling cream DAB 6 (stabili...	15066601	2
Caelo	Cooling cream DAB	1506601	1
Caelo	Cooling cream DAB	15066601	13
Sanacorp	Cooling cream DAB 6 (stabili...	15066601	2
Alliance	Cooling cream DAB	15066601	1
Caelo	Cooling cream DAB	15066701	4
Noweda	Cooling cream DAB	15066601	1
Caelo	Cooling cream DAB	15066801	5
Caelo	Cooling cream DAB	15066801test	1
Caelo	Cooling cream DAB	15066901	1
Caelo	Cooling cream DAB	15069201	1
cealo	Cooling cream DAB	15066901	1
Caelo	Cooling cream DAB	1513101	1
Caelo	Cooling cream DAB	1509201	1
Caelo	Cooling cream DAB	15152901	4
Caelo	Cooling cream DAB 6 (stabili...	15153001	1
Caelo	Cooling cream DAB	15153001	3
Caelo	Cooling cream DAB 6 (stabili...	15153101	1
Caelo	Cooling cream DAB	15153201	7
Caelo	Cooling cream DAB 6 (stabili...	15153201	1
Caelo	Cooling cream DAB	15153101	3
Caelo	Cooling cream DAB	15153202	1
Caelo	Cooling cream DAB 6 (stabili...	15153301	2
Caelo	Cooling cream DAB	15153301	3
Caelo	Cooling cream DAB	15153401	1
Caelo	Cooling cream DAB	15153501	17
Caelo	Cooling cream DAB 6 (stabili...	15153501	2
Caelo	Cooling cream DAB 6 (stabili...	15228401	2
Caelo	Cooling cream DAB	15228401	4
Cealo/ Noweda	Cooling cream DAB	15228501	1
Caesar&Loretz	Cooling cream DAB	15228501	2
Caelo	Cooling cream DAB 6 (stabili...	15228502	3
Caelo	Cooling cream DAB	15228502	1
Caelo	Cooling cream DAB 6 (stabili...	15228601	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Cae/Phönix	Cooling cream DAB	15228601	1
sanacorp	Cooling cream DAB	15228601	1
Rondell Apotheke	Cooling cream DAB	15228601	1
Caesar & Loretz GmbH	Cooling cream DAB	15228601	1
Caelo	Cooling cream DAB	15228701	8
Chiesi	Cooling cream DAB	15228701	1
Caelo	Cooling cream DAB	15228801	8
Caelo	Cooling cream DAB 6 (stabili...	152327901	1
Caelo	Cooling cream DAB	1523609	1
Caelo	Cooling cream DAB	15236701	1
Fagron	Cooling cream DAB	15236801	1
Caelo	Cooling cream DAB	15236801	1
Caelo	Cooling cream DAB 6 (stabili...	15236801	1
Caelo	Cooling cream DAB	15236901	19
Caelo	Cooling cream DAB 6 (stabili...	15236901	4
Hubertus Apotheke	Cooling cream DAB 6 (stabili...	15236901	1
Caelo	Cooling cream DAB	15237101	7
Caelo	Cooling cream DAB	15327701	6
Caelo	Cooling cream DAB	15327801	3
Caelo	Cooling cream DAB	15327901	3
Caelo	Cooling cream DAB 6 (stabili...	15327901	1
Caelo	Cooling cream DAB	15328101	2
Caelo	Cooling cream DAB	15413501	3
Caelo	Cooling cream DAB	15413601	4
Caelo	Cooling cream DAB 6 (stabili...	15413501	1
Caelo	Cooling cream DAB 6 (stabili...	15413601	2
Caelo	Cooling cream DAB	15413701	4
Sanacorp	Cooling cream DAB	15413901	1
Caelo	Cooling cream DAB	15417201	6
Caelo	Cooling cream DAB 6 (stabili...	15413901	1
Noweda	Cooling cream DAB	15417201	1
Caelo	Cooling cream DAB	15417301	2
Caelo	Cooling cream DAB	15417601	4
Caelo	Cooling cream DAB 6 (stabili...	15417401	1
Caelo	Cooling cream DAB 6 (stabili...	15417501	1
Caesar & Loretz GmbH	Cooling cream DAB	15434401	1
Cordes/Phönix	Cooling cream DAB	15434402	1
Caelo	Cooling cream DAB 6 (stabili...	15434402	1
Cealo/AHD	Cooling cream DAB	15434402	1
phönix	Cooling cream DAB	15435101	1
Caelo	Cooling cream DAB 6 (stabili...	1547501	1
Caelo	Cooling cream DAB	1547501	1
Caelo	Cooling cream DAB	16071505	1
Cealo/AHD	Cooling cream DAB 6 (stabili...	16011901	1
Caelo	Cooling cream DAB 6 (stabili...	1608601	1
Caelo	Cooling cream DAB 6 (stabili...	16123501	1
gehe	Cooling cream DAB	17189701	1
phönix	Cooling cream DAB	17033001	1
Bombastus	Cooling cream DAB		1
Bombastus	Cooling cream DAB	19032014A	1
Caelo	Cooling cream DAB	21071401	1
Caelo	Cooling cream DAB	19051506	1
Caelo	Cooling cream DAB	23121412	1
Caelo	Cooling cream DAB	19051505	1
Bombastus	Cooling cream DAB	2405212541	1
Caelo	Cooling cream DAB	27061201	1
Bombastus	Cooling cream DAB	271726	1
Bombastus	Cooling cream DAB	282021	2
Bombastus	Cooling cream DAB	284264	1
Caelo	Cooling cream DAB	283118	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Cooling cream DAB	283119	1
Bombastus	Cooling cream DAB	284299	2
Bombastus	Cooling cream DAB	284298	1
Bombastus	Cooling cream DAB	284302	3
Bombastus	Cooling cream DAB	284303	3
Bombastus	Cooling cream DAB	285063	1
Bombastus	Cooling cream DAB	285068	2
Bombastus	Cooling cream DAB	285069	3
Bombastus	Cooling cream DAB	285065	1
Bombastus	Cooling cream DAB	285498	1
gehe	Cooling cream DAB	285498	1
Bombastus	Cooling cream DAB	285634	2
Bombastus	Cooling cream DAB	285502	2
Bombastus	Cooling cream DAB	285504	2
Bombastus	Cooling cream DAB	285503	1
Bombastus	Cooling cream DAB	285635	2
Bombastus	Cooling cream DAB	285636	1
Bombastus	Cooling cream DAB	285637	1
Bombastus	Cooling cream DAB	285638	2
Bombastus	Cooling cream DAB	285639	1
Bombastus	Cooling cream DAB	285644	2
Bombastus	Cooling cream DAB	285645	2
Bombastus	Cooling cream DAB	286562	1
Bombastus	Cooling cream DAB	286564	3
Gehe	Cooling cream DAB	286565	1
Fontapharm/Phönix	Cooling cream DAB	286565	1
Bombastus	Cooling cream DAB	286648	3
Bombastus	Cooling cream DAB	286567	1
Bombastus	Cooling cream DAB	286654	4
Bombastus	Cooling cream DAB	286650	3
Bombastus	Cooling cream DAB	286653	2
Bombastus	Cooling cream DAB	287485	2
Bombastus	Cooling cream DAB	287483	1
Phönix	Cooling cream DAB	287486	1
Bombastus	Cooling cream DAB	287542	1
Bombastus	Cooling cream DAB	287544	2
Bombastus	Cooling cream DAB	287539	1
Caelo	Cooling cream DAB	287544	1
Bombastus	Cooling cream DAB	287543	1
Bombastus	Cooling cream DAB	288107	1
Bombastus	Cooling cream DAB	288109	2
Bombastus	Cooling cream DAB	288112	2
Bombastus	Cooling cream DAB	288110	1
Bombastus	Cooling cream DAB	288152	3
Bombastus	Cooling cream DAB	288158	1
Gehe	Cooling cream DAB	289187	1
Bombastus	Cooling cream DAB	288159	1
Bombastus	Cooling cream DAB	289188	2
Bombastus	Cooling cream DAB	289187	1
Bombastus	Cooling cream DAB	289189	1
Phönix	Cooling cream DAB	289188	1
Bombastus	Cooling cream DAB	289251	1
Bombastus	Cooling cream DAB	289253	2
Bombastus	Cooling cream DAB	289191	3
Bombastus	Cooling cream DAB	289256	1
Bombastus	Cooling cream DAB	289258	1
Bombastus	Cooling cream DAB	289260	3
Caelo	Cooling cream DAB	289265	1
Bombastus	Cooling cream DAB	290380	2
Bombastus	Cooling cream DAB	290383	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Cooling cream DAB	290396	1
Bombastus	Cooling cream DAB	290483	2
Gehe	Cooling cream DAB	290485	1
Bombastus	Cooling cream DAB 6 (stabili...	290486	1
Beiersdorf / E	Cooling cream DAB	290486	1
Bombastus	Cooling cream DAB	290487	1
Bombastus	Cooling cream DAB	290489	1
Bombastus	Cooling cream DAB	290490	1
Caelo	Cooling cream DAB	29081509	1
Bombastus	Cooling cream DAB	29072015A	1
Bombastus	Cooling cream DAB	291558	1
Bombastus	Cooling cream DAB	291552	1
Bombastus	Cooling cream DAB	291555	1
Bombastus	Cooling cream DAB	291618	1
Bombastus	Cooling cream DAB	291620	3
Bombastus	Cooling cream DAB	292103	1
Bombastus	Cooling cream DAB	292106	3
Bombastus	Cooling cream DAB	292202	1
Bombastus	Cooling cream DAB	292107	1
Bombastus	Cooling cream DAB	292280	1
Bombastus	Cooling cream DAB	292102	1
Bombastus	Cooling cream DAB	292281	1
Bombastus	Cooling cream DAB	292930	2
Bombastus	Cooling cream DAB	292936	1
Bombastus	Cooling cream DAB	292997	4
Bombastus	Cooling cream DAB	292994	1
Bombastus	Cooling cream DAB	293616	2
Bombastus	Cooling cream DAB	293622	2
Bombastus	Cooling cream DAB	293621	3
Gehe	Cooling cream DAB	293677	1
Bombastus	Cooling cream DAB	293747	1
Caelo	Cooling cream DAB	293748	1
Bombastus	Cooling cream DAB	293749	1
Bombastus	Cooling cream DAB	293750	1
Bombastus	Cooling cream DAB	295019	2
Bombastus	Cooling cream DAB	295021	1
Bombastus	Cooling cream DAB	295120	1
Bombastus	Cooling cream DAB	295183	2
EB-Ho	Cooling cream DAB	295184	1
Bombastus	Cooling cream DAB	295189	1
Caelo	Cooling cream DAB	295123	1
Bombastus	Cooling cream DAB	295184	2
Bombastus	Cooling cream DAB	295185	1
Bombastus	Cooling cream DAB	295190	1
Bombastus	Cooling cream DAB	295192	1
Bombastus	Cooling cream DAB 6 (stabili...	295191	1
Bombastus	Cooling cream DAB	295760	2
Bombastus	Cooling cream DAB	295756	1
Bombastus	Cooling cream DAB	295761	1
Caelo	Cooling cream DAB	295762	1
Bombastus	Cooling cream DAB	296312	3
Bombastus	Cooling cream DAB	296316	1
bomb	Cooling cream DAB	296316	1
Bombastus	Cooling cream DAB	296317	1
Bombastus	Cooling cream DAB	297097	2
Bombastus	Cooling cream DAB	297099	1
Bombastus	Cooling cream DAB	297102	1
Bombastus	Cooling cream DAB	297205	2
Bombastus	Cooling cream DAB	297206	1
Bombastus	Cooling cream DAB	297209	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Cooling cream DAB	297210	2
Bombastus	Cooling cream DAB	297211	1
Bombastus	Cooling cream DAB	297214	2
Bombastus	Cooling cream DAB	298118	1
Bombastus	Cooling cream DAB	298121	1
Bombastus	Cooling cream DAB 6 (stabili...	298319	1
Bombastus	Cooling cream DAB	298879	1
Bombastus	Cooling cream DAB	298881	3
apomix Fiebig	Cooling cream DAB	299362	1
Bombastus	Cooling cream DAB	299361	1
Bombastus	Cooling cream DAB	299973	1
Bombastus	Cooling cream DAB	299976	2
Caelo	Cooling cream DAB	30011505	1
Caelo	Cooling cream DAB	3012M-02541	1
Caelo	Cooling cream DAB	30091402	1
Bombastus	Cooling cream DAB	300724	1
Henry Lamotte/GEHE	Cooling cream DAB	5163700	1
PHOENIX am 06.06.2012	Cooling cream DAB	5163800	1
Caelo	Cooling cream DAB 6 (stabili...	5328001	1
Caelo	Cooling cream DAB	5600300	1
Lamotte/ Phönix	Cooling cream DAB	5600300	1
Henry Lamotte Oils/GEHE	Cooling cream DAB	5600300	1
Lamotte Sanacorp	Cooling cream DAB	5496200	3
Henry Lamotte Olis GmbH...	Cooling cream DAB	5600300	2
Lamotte Sanacorp	Cooling cream DAB	5600300	1
Sanacorp/Henry Lamotte	Cooling cream DAB	6000083005	1
Henry Lamotte Olis GmbH	Cooling cream DAB	5735600	1
Henry Lamotte Oils/GEHE	Cooling cream DAB	5691200	1
Caelo	Cooling cream DAB	5735600	1
Henry Lamotte Oils/GEHE	Cooling cream DAB	5735600	1
Lamotte/AHD	Cooling cream DAB	6001092001	1
Anzag	Cooling cream DAB	6002062002	1
Bombastuts/Gehe	Cooling cream DAB	60082802	1
Henry Lamotte Oils	Cooling cream DAB	6003682010	1
Sanacorp	Cooling cream DAB	77/AKD10	1
Bombastus	Cooling cream DAB	91709	1
Sanacorp	Cooling cream DAB	77/AKD6	1
Fagron	Cooling cream DAB	B15A30T07018676	1
Phoenix	Cooling cream DAB	269094	1
Gehe	Cooling cream DAB	270314	1
Bombastus	Cooling cream DAB	271541	1
Bombastus	Cooling cream DAB	270315	1
Bombastus	Cooling cream DAB	269092	1
Phönix	Cooling cream DAB	271541	1
Bombastus	Cooling cream DAB	271623	1
Bombastus	Cooling cream DAB	271624	1
Phoenix 12.11.12	Cooling cream DAB	271623	1
Caelo	Cooling cream DAB	272479	2
Bombastus	Cooling cream DAB	271627	1
Bombastus	Cooling cream DAB	272353	1
Caelo	Cooling cream DAB	271631	1
Bombastus	Cooling cream DAB	272349	1
Bombastus	Cooling cream DAB	272480	1
Bombastus	Cooling cream DAB	272482	1
Anzag 25.10.12	Cooling cream DAB	272484	1
Bombastus	Cooling cream DAB	272485	1
Anzag, 25.10.2012	Cooling cream DAB	272486	1
Bombastus	Cooling cream DAB	272564	2
Bombastus	Cooling cream DAB	272573	2
Bombastus	Cooling cream DAB	272572	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	273441	3
Bombastus	Cooling cream DAB	273443	1
Bombastus	Cooling cream DAB	273444	1
Bombastus	Cooling cream DAB	273503	1
Bombastus	Cooling cream DAB	273499	1
Noweda	Cooling cream DAB	274043	1
Bombastus	Cooling cream DAB	273501	1
ANZAG	Cooling cream DAB	273504	1
Bombastus	Cooling cream DAB	274044	2
Gehe	Cooling cream DAB	274044	2
Bombastus	Cooling cream DAB	274077	1
Sanacorp	Cooling cream DAB	274077	1
Bombastus	Cooling cream DAB	274079	1
Bombastus	Cooling cream DAB	274080	2
Bombastus	Cooling cream DAB	274083	2
Caelo	Cooling cream DAB	6061303	1
Caelo	Cooling cream DAB	2081201	1
Bombastus	Cooling cream DAB	513A-02541	1
Caelo	Cooling cream DAB	6111201	1
Caelo	Cooling cream DAB	613I-02541	1
Caelo	Cooling cream DAB	6111202	1
Bombastus	Cooling cream DAB	614U2541	1
Caelo	Cooling cream DAB	7021301	1
Caelo	Cooling cream DAB	7091201	1
Caelo	Cooling cream DAB	913I-02541	1
Caelo	Cooling cream DAB	8041304	1
Caelo	Cooling cream DAB	10081209	1
Caelo	Cooling cream DAB	10091306	1
Caelo	Cooling cream DAB	8111201	1
Bombastus	Cooling cream DAB	1109212541	1
Caelo	Cooling cream DAB	10081210	1
Caelo	Cooling cream DAB	11269101	1
Caelo	Cooling cream DAB	1113Q-02541	1
Caelo/Gehe	Cooling cream DAB	12052002	2
Caelo	Cooling cream DAB	12052001	1
Caelo	Cooling cream DAB	120522/1	1
Bombastus	Cooling cream DAB	12082013A	1
Caelo	Cooling cream DAB	12089201	1
Caelo	Cooling cream DAB	12089301	2
Caelo	Cooling cream DAB	12089401	1
Caelo	Cooling cream DAB	12089402	2
Phönix Gotha	Cooling cream DAB	12089501	1
Caelo	Cooling cream DAB	12090201	3
Caelo	Cooling cream DAB	12090501	2
VDL	Cooling cream DAB	12090501	1
Phoenix, 04.09.2012	Cooling cream DAB	12090601	2
Caelo	Cooling cream DAB	1213M-02541	2
Caelo	Cooling cream DAB	12090601	1
Noweda	Cooling cream DAB	12174301	3
Caelo	Cooling cream DAB	12174301	2
Caelo	Cooling cream DAB	12174401	4
Caesar & Loretz GmbH, 4. . .	Cooling cream DAB	12174401	2
Caelo	Cooling cream DAB	12174601	2
Caelo	Cooling cream DAB	12174501	1
Caelo	Cooling cream DAB	12178801	2
Phönix, 30.10.2012	Cooling cream DAB	12178801	1
Caesar & Loretz GmbH	Cooling cream DAB	12178901	1
Caelo	Cooling cream DAB	12178902	3
Caelo	Cooling cream DAB	12178901	1
Großhandel	Cooling cream DAB	12178902	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	12179001	3
Erika apotheke	Cooling cream DAB	12182701	1
Caelo	Cooling cream DAB	12182701	5
Caesar & Loretz GmbH, 4...	Cooling cream DAB	12182801	2
Caelo	Cooling cream DAB	12182801	3
Caelo	Cooling cream DAB	12187201	3
Sanacorp	Cooling cream DAB	12187201	1
Caelo	Cooling cream DAB	12187301	4
Kehr	Cooling cream DAB	12187301	1
Sanacorp	Cooling cream DAB	12187301	1
Phönix 26.10.12	Cooling cream DAB	12187301	1
Caelo	Cooling cream DAB	12253001	3
Caesar & Loretz GmbH D-...	Cooling cream DAB	12253001	1
Caelo	Cooling cream DAB	12253101	4
Caelo	Cooling cream DAB	12253201	3
Großhandel	Cooling cream DAB	12253201	1
Caelo	Cooling cream DAB	12253301	2
Caelo	Cooling cream DAB	12253401	2
Caelo	Cooling cream DAB	12290102	1
Caelo	Cooling cream DAB	12291301	2
Caelo	Cooling cream DAB	12291901	2
Anzag 21.02.2013	Cooling cream DAB	12291901	1
Caelo	Cooling cream DAB	12292102	17
Caelo	Cooling cream DAB	12292201	6
Großhandel	Cooling cream DAB	12292201	1
phönix 22.01.2013	Cooling cream DAB	12292201	2
Caelo	Cooling cream DAB	12292301	5
phönix	Cooling cream DAB	12292201	1
Caesar & Loretz GmbH, 4...	Cooling cream DAB	12292301	1
Henry Lamotte Oils GmbH	Cooling cream DAB	12292301	3
Noweda	Cooling cream DAB	12292301	1
Gehe	Cooling cream DAB	12292401	1
Caesar & Loretz GmbH	Cooling cream DAB	12306001	1
Kehr	Cooling cream DAB	12312002	1
Caelo	Cooling cream DAB	12312101	2
phönix	Cooling cream DAB	12312102	1
Caelo	Cooling cream DAB	12312201	5
Caelo	Cooling cream DAB	12312301	1
Caelo	Cooling cream DAB	12312401	9
Sanacorp	Cooling cream DAB	12312401	1
Caelo	Cooling cream DAB	12312501	3
Caelo	Cooling cream DAB	12312601	3
Caelo	Cooling cream DAB	12312701	3
Caelo	Cooling cream DAB	12312801	3
Großhandel	Cooling cream DAB	12312801	1
Gehe	Cooling cream DAB	12312801	1
Sanacorp	Cooling cream DAB	12312801	1
Caelo	Cooling cream DAB	12364701	2
Caelo	Cooling cream DAB	12364801	4
GEHE	Cooling cream DAB	12364801	1
Caelo	Cooling cream DAB	12364802	1
Caelo	Cooling cream DAB	1287201	1
Caelo	Cooling cream DAB	1289501	1
Caelo	Cooling cream DAB	1290301	1
Fagron	Cooling cream DAB	12G25-T02	2
Noweda	Cooling cream DAB	12992301	1
Caelo	Cooling cream DAB	13004101	5
Caelo	Cooling cream DAB	13004201	8
cael	Cooling cream DAB	13004201	1
Caelo	Cooling cream DAB	13004401	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	13004402	6
Caelo	Cooling cream DAB	13004301	3
Caesar/Lorenz	Cooling cream DAB	13004401	1
Kehr	Cooling cream DAB	13004501	1
Caelo	Cooling cream DAB	13004501	7
Großhandel	Cooling cream DAB	13004601	1
Caelo	Cooling cream DAB	13004601	4
Gehe	Cooling cream DAB	13004701	1
calo	Cooling cream DAB	13004701	2
Caelo	Cooling cream DAB	13026201	1
Bombastus	Cooling cream DAB	13012014A	1
Caelo	Cooling cream DAB	13051401	5
Caelo	Cooling cream DAB	13047201	1
Noweda	Cooling cream DAB	13051501	1
Caelo	Cooling cream DAB	13051501	9
Caelo	Cooling cream DAB	13051602	4
Phönix	Cooling cream DAB	13051602	1
Caelo	Cooling cream DAB	13059501	2
Phönix	Cooling cream DAB	13059501	1
Großhandel	Cooling cream DAB	13059501	1
Caelo	Cooling cream DAB	13059601	10
Sanacorp	Cooling cream DAB	13059601	1
Caesar & Loretz GmbH, 4...	Cooling cream DAB	13059601	1
Großhandel	Cooling cream DAB	13059601	1
Caelo	Cooling cream DAB	13059701	7
Phönix	Cooling cream DAB	13059901	1
Caelo	Cooling cream DAB	13059902	3
Caelo	Cooling cream DAB	13068101	3
Caelo	Cooling cream DAB	13121001	10
Caelo	Cooling cream DAB	13121002	4
Lamotte/Phönix	Cooling cream DAB	1312101	3
Caelo	Cooling cream DAB	13121201	7
Caelo	Cooling cream DAB	13121301	5
Caesar & Loretz GmbH, 4...	Cooling cream DAB	13121401	2
Caelo	Cooling cream DAB	13121401	10
Caelo	Cooling cream DAB	13121501	1
Caelo	Cooling cream DAB	13121503	1
Caelo	Cooling cream DAB	13121601	2
Caelo	Cooling cream DAB	13121701	7
Caesar & Loretz GmbH, 4...	Cooling cream DAB	13121701	2
Großhandel	Cooling cream DAB	13121701	1
Caelo	Cooling cream DAB	1312601	1
Caelo	Cooling cream DAB	13174901	1
Caelo	Cooling cream DAB	13175001	2
Anzag	Cooling cream DAB	13175101	1
Caesar & Loretz GmbH, 4...	Cooling cream DAB	13175101	4
Caelo	Cooling cream DAB	13191301	2
Caelo	Cooling cream DAB	13191401	11
Gehe	Cooling cream DAB	13191501	1
Caelo	Cooling cream DAB	13191601	6
Großhandel	Cooling cream DAB	13191501	1
Großhandel	Cooling cream DAB	13191801	1
Maros/Phoenix	Cooling cream DAB	13191901	1
Caelo	Cooling cream DAB	13192001	7
Caesar & Loretz GmbH, 4...	Cooling cream DAB	13192001	1
Phönix	Cooling cream DAB	13192101	1
Sanacorp	Cooling cream DAB	13192101	1
Caelo	Cooling cream DAB	13192201	1
Caelo	Cooling cream DAB	13225701	2
Caelo	Cooling cream DAB	1322601	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	13225801	1
AHD, 08.11.13	Cooling cream DAB	13225801	1
Caelo	Cooling cream DAB	13207201	1
Caelo	Cooling cream DAB	13226301	8
Caelo	Cooling cream DAB	13226302	6
Phönix	Cooling cream DAB	13226302	1
Sanacorp	Cooling cream DAB	13226501	1
Lamotte	Cooling cream DAB	13226601	1
Caelo	Cooling cream DAB	13226602	1
Caelo	Cooling cream DAB	13226701	4
Sanacorp	Cooling cream DAB	13226701	1
Klenk	Cooling cream DAB	13226801	1
Caelo	Cooling cream DAB	13226901	8
Caelo	Cooling cream DAB	13227001	3
Gehe am 03.02.2014	Cooling cream DAB	1322701	1
Großhandel	Cooling cream DAB	13227001	3
Caelo	Cooling cream DAB	13227101	6
Caelo	Cooling cream DAB	13227201	8
Euro OTC	Cooling cream DAB	13227101	1
Caelo	Cooling cream DAB	13227301	4
Caelo	Cooling cream DAB	13227401	9
Bombastus	Cooling cream DAB	13315101	2
Caelo	Cooling cream DAB	13315102	3
Caesar & Loretz GmbH	Cooling cream DAB	13315101	1
gehe	Cooling cream DAB	13315102	1
Noweda	Cooling cream DAB	13315101	1
Gehe	Cooling cream DAB	13315201	1
Caesar & Loretz GmbH	Cooling cream DAB	13315201	1
Caelo	Cooling cream DAB	13315401	1
Caesar & Loretz GmbH, 4...	Cooling cream DAB	13315601	1
Bombastus	Cooling cream DAB	13315701	1
Caelo	Cooling cream DAB	13315901	1
Caelo	Cooling cream DAB	133227401	1
Caelo	Cooling cream DAB	13327101	1
Sanacorp	Cooling cream DAB	13401801	1
Caelo	Cooling cream DAB	13402001	5
Sanacorp	Cooling cream DAB	13402001	1
Fagron	Cooling cream DAB	13D19-T01-002941	1
Fagron	Cooling cream DAB	13E10-T01-003539	2
Caelo	Cooling cream DAB	14055703	2
Caelo	Cooling cream DAB	1408601	1
Caelo	Cooling cream DAB	15011304	1
Bombastus	Cooling cream DAB	151113a	1
Caelo	Cooling cream DAB	1543512	1
Fagron	Cooling cream DAB	15A30-T07-018974	1
Caelo	Cooling cream DAB	16081303	1
Caelo	Cooling cream DAB	16101201	1
Bombastus	Cooling cream DAB	2112212541	1
Caelo	Cooling cream DAB	17121301	1
Caelo	Cooling cream DAB	18121205	1
Caelo	Cooling cream DAB	22021305	1
Caelo	Cooling cream DAB	25021414	1
Bombastus	Cooling cream DAB	2413Q02541	1
Caelo	Cooling cream DAB	22051307	1
Bombastus	Cooling cream DAB	22487	1
Caelo	Cooling cream DAB	22061301	1
Caelo	Cooling cream DAB	25091308	1
Caelo	Cooling cream DAB	26071307	1
Caelo	Cooling cream DAB	27111201	1
Bombastus	Cooling cream DAB	2713Q2541	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Cooling cream DAB	272571	1
Gehe	Cooling cream DAB	274084	1
Bombastus	Cooling cream DAB	274042	1
Bombastus	Cooling cream DAB	274785	1
Bombastus	Cooling cream DAB	274791	2
Gehe	Cooling cream DAB	274909	2
Bombastus	Cooling cream DAB	274914	1
Anzag	Cooling cream DAB	274911	1
Bombastus	Cooling cream DAB	275663	1
Bombastus	Cooling cream DAB	275664	2
Caelo	Cooling cream DAB	275666	1
Bombastus	Cooling cream DAB	275670	1
Bombastus	Cooling cream DAB	275668	1
Bombastus	Cooling cream DAB	276282	1
Bombastus	Cooling cream DAB	276281	1
Gehe	Cooling cream DAB	276283	3
Bombastus	Cooling cream DAB	276288	1
Noweda	Cooling cream DAB	276284	1
Bombastus	Cooling cream DAB	276418	4
Bombastus	Cooling cream DAB	276284	1
Bombastus	Cooling cream DAB	276422	2
Bombastus	Cooling cream DAB	276419	2
Bombastus	Cooling cream DAB	277275	2
Bombastus	Cooling cream DAB	277363	1
Alliance Healthcare	Cooling cream DAB	277366	1
Bombastus	Cooling cream DAB	277525	2
Alliance Healthcare	Cooling cream DAB	277368	1
Bombastus	Cooling cream DAB	277530	2
Bombastus	Cooling cream DAB	277526	1
Bombastus	Cooling cream DAB	277371	1
Bombastus	Cooling cream DAB	277533	1
Fagron	Cooling cream DAB	277533	1
Alliance Healthcare	Cooling cream DAB	277534	1
Bombastus	Cooling cream DAB	277534	2
Bombastus	Cooling cream DAB	278337	1
Bombastus	Cooling cream DAB	278338	1
Bombastus	Cooling cream DAB	278439	2
Bombastus	Cooling cream DAB	278437	2
Bombastus	Cooling cream DAB	278436	1
Bombastus	Cooling cream DAB	278440	3
Bombastus	Cooling cream DAB	278342	1
Bombastus	Cooling cream DAB	278441	1
Alliance Healthcare	Cooling cream DAB	278444	1
Phönix	Cooling cream DAB	278516	2
Gehe	Cooling cream DAB	278518	1
Bombastus	Cooling cream DAB	278525	3
Alliance Healthcare	Cooling cream DAB	278523	1
Bombastus	Cooling cream DAB	278674	1
Bombastus	Cooling cream DAB	279680	2
Gehe	Cooling cream DAB	279681	2
Bombastus	Cooling cream DAB	279681	1
Bombastus	Cooling cream DAB	279682	2
Bombastus	Cooling cream DAB	279750	2
Bombastus	Cooling cream DAB	279751	1
Bombastus	Cooling cream DAB	279755	3
Bombastus	Cooling cream DAB	279765	2
Bombastus	Cooling cream DAB	279766	1
Bombastus	Cooling cream DAB	281033	1
Bombastus	Cooling cream DAB	281030	1
Gehe	Cooling cream DAB	281036	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cooling cream DAB	28101302	1
Caelo	Cooling cream DAB	28021306	1
Bombastus	Cooling cream DAB	281138	1
Bombastus	Cooling cream DAB	281134	1
Sanacorp	Cooling cream DAB	281037	1
Bombastus	Cooling cream DAB	281139	3
Bombastus	Cooling cream DAB	281140	1
Caelo	Cooling cream DAB	28121203	1
Bombastus	Cooling cream DAB	281820	1
Bombastus	Cooling cream DAB	281822	1
Bombastus	Cooling cream DAB	281824	1
Bombastus	Cooling cream DAB	281912	1
Bombastus	Cooling cream DAB	282017	1
Bombastus	Cooling cream DAB	281919	1
Bombastus	Cooling cream DAB	282019	2
Caesar & Loretz GmbH	Cooling cream DAB	282019	1
Gehe	Cooling cream DAB	282020	1
Caelo	Cooling cream DAB	282020	1
Gehe	Cooling cream DAB	283053	1
Gehe	Cooling cream DAB	282022	2
Bombastus	Cooling cream DAB	283057	1
Bombastus	Cooling cream DAB	283051	1
Bombastus	Cooling cream DAB	283054	1
Bombastus	Cooling cream DAB	283058	1
Gehe	Cooling cream DAB	283114	1
Bombastus	Cooling cream DAB	283116	1
Bombastus	Cooling cream DAB	283117	1
Bombastus	Cooling cream DAB	283118	1
Bombastus	Cooling cream DAB	284263	1
Gehe	Cooling cream DAB	284300	1
Bombastus	Cooling cream DAB	284300	3
Gehe	Cooling cream DAB	285066	1
Bombastus	Cooling cream DAB	285070	1
Euro OTC	Cooling cream DAB	286563	1
Bombastus	Cooling cream DAB	286652	1
Bombastus	Cooling cream DAB	297213	1
Caelo	Cooling cream DAB	30081302	1
Caelo	Cooling cream DAB	31011301	1
Caelo	Cooling cream DAB	3963/V12/14	1
Caelo	Cooling cream DAB	3413E-02541	1
Caelo	Cooling cream DAB	3964/V01/15	1
Bombastus	Cooling cream DAB	4013Q2541	1
Caelo	Cooling cream DAB	4313I-02541	1
Anzag/	Cooling cream DAB	5316800	1
Henry Lamotte/GEHE	Cooling cream DAB	5316800	1
Caelo	Cooling cream DAB	5132	1
phönix 17.1.13	Cooling cream DAB	5316900	1
Noweda/Lamotte Oils	Cooling cream DAB	5321500	1
Henry Lamotte	Cooling cream DAB	5321500	1
Henry Lamotte/GEHE	Cooling cream DAB	5321500	1
Caelo	Cooling cream DAB	5321500	1
Henry Lamotte	Cooling cream DAB	5407300	2
Henry Lamotte Oils/GEHE	Cooling cream DAB	5407600	1
Anzag/Lamotte	Cooling cream DAB	5407300	1
Noweda/Lamotte	Cooling cream DAB	5496200	2
Henry Lamotte	Cooling cream DAB	5496300	1
Bombastus	Cooling cream DAB	7022014B	1
Bombastus	Cooling cream DAB	60082802	1
Caelo	Cooling cream DAB	72/AD2	1
Bombastus	Cooling cream DAB	77309	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Cooling cream DAB	caelo	1
Phönix Gotha	Cooling cream DAB	SRTM1104	1
Bombastus	Cooling cream DAB	92086	1
Caelo	Cooling cream DAB 6 (stabili...	12179001	1
Caelo	Cooling cream DAB 6 (stabili...	213E-02541	1
Caelo	Cooling cream DAB 6 (stabili...	12182701	2
Gehe	Cooling cream DAB 6 (stabili...	12312601	1
Caelo	Cooling cream DAB 6 (stabili...	12312701	1
Noweda	Cooling cream DAB 6 (stabili...	13026301	2
Caelo	Cooling cream DAB 6 (stabili...	12364802	1
Caelo	Cooling cream DAB 6 (stabili...	13047201	1
Caelo	Cooling cream DAB 6 (stabili...	13051501	2
Caelo	Cooling cream DAB 6 (stabili...	13051401	1
Caelo	Cooling cream DAB 6 (stabili...	13121401	1
Caelo	Cooling cream DAB 6 (stabili...	13120901	3
Caelo	Cooling cream DAB 6 (stabili...	13175101	1
Caelo	Cooling cream DAB 6 (stabili...	13191901	1
Caelo	Cooling cream DAB 6 (stabili...	13175301	3
Caelo	Cooling cream DAB 6 (stabili...	13192101	5
Caelo	Cooling cream DAB 6 (stabili...	14086301	1
Gehe	Cooling cream DAB 6 (stabili...	12003002	1
Caelo	Cooling cream DAB 6 (stabili...	12206101	1
Caelo	Cooling cream DAB 6 (stabili...	20130912-3	1
Gehe	Cooling cream DAB 6 (stabili...	273502	1
anzag, 21.01.2013	Cooling cream DAB 6 (stabili...	12312301	1
Caelo	Cooling cream DAB 6 (stabili...	12312501	1
Caelo	Cooling cream DAB 6 (stabili...	12312401	1
Caelo	Cooling cream DAB 6 (stabili...	12364801	1
Caelo	Cooling cream DAB 6 (stabili...	13004401	1
Caelo	Cooling cream DAB 6 (stabili...	1302601	1
Caelo	Cooling cream DAB 6 (stabili...	13051601	2
Caelo	Cooling cream DAB 6 (stabili...	13092901	1
Caelo	Cooling cream DAB 6 (stabili...	13121201	5
Sanacorp	Cooling cream DAB 6 (stabili...	13121601	1
Caelo	Cooling cream DAB 6 (stabili...	13121701	1
Caelo	Cooling cream DAB 6 (stabili...	13174901	1
Caelo	Cooling cream DAB 6 (stabili...	13191501	1
Sanacorp	Cooling cream DAB 6 (stabili...	13191301	1
gehe 14.08.2013	Cooling cream DAB 6 (stabili...	13191401	1
Caelo	Cooling cream DAB 6 (stabili...	13191701	1
Caelo	Cooling cream DAB 6 (stabili...	13192001	1
Caelo	Cooling cream DAB 6 (stabili...	13226302	1
Caelo	Cooling cream DAB 6 (stabili...	13226601	2
Caelo	Cooling cream DAB 6 (stabili...	13226801	1
Caelo	Cooling cream DAB 6 (stabili...	14129801	1
Caelo	Cooling cream DAB 6 (stabili...	14085901	1
Caelo	Cooling cream DAB 6 (stabili...	14184601	1
Caelo	Cooling cream DAB	14268002	1
	Cooling cream DAB	14328801	1
Caelo	Cooling cream DAB	16000401	1
sanacorp, 9,51EUR, 27.09...	Cooling cream DAB	16172001	1
Caelo	Cooling cream DAB 6 (stabili...	17000601	1
Bombastus	Cooling cream DAB	285499	2
Bombastus	Cooling cream DAB	289265	1
Bombastus	Cooling cream DAB	291619	1
Bombastus	Cooling cream DAB	295124	1
Lamotte/Spangro	Cooling cream DAB	5735700	1
Lamotte/Spangro	Cooling cream DAB	5861100	12
Lamotte/Spangro	Cooling cream DAB	6000744001	3
Lamotte/Spangro	Cooling cream DAB	6001092002	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte / AHD	Cooling cream DAB	8003478004	1
Gutenberg-Apotheke	Cooling cream DAB	84187	1
Caelo	Cooling cream DAB	P158/14	1
Bombastus	Cooling cream DAB	270310	1
Anzag	Cooling cream DAB	272348	1
Caelo	Cooling cream DAB	11323501	1
Caelo	Cooling cream DAB	12312001	1
Sanacorp	Cooling cream DAB	13004401	2
Caelo	Cooling cream DAB	1312091	1
Sanacorp	Cooling cream DAB	13121001	1
A.H.D.	Cooling cream DAB	13226301	1
Caleo/Phönix	Cooling cream DAB	14086501	1
Anzag	Cooling cream DAB	274042	1
Bombastus	Cooling cream DAB	274783	1
Anzag	Cooling cream DAB	272347	1
Bombastus	Cooling cream DAB	275665	1
Bombastus	Cooling cream DAB	277367	1
Bomba/AHD	Cooling cream DAB	277269	1
Bombastus	Cooling cream DAB	275667	1
Bombastus	Cooling cream DAB	277369	1
Bombastus	Cooling cream DAB	278335	3
Bombastus	Cooling cream DAB	278341	2
Bombastus	Cooling cream DAB	278438	1
Bombastus	Cooling cream DAB	279676	2
Bombastus	Cooling cream DAB	278442	1
Pharmachem/AHD	Cooling cream DAB	279680	1
Bomba/AHD	Cooling cream DAB	281031	2
Caelo	Cooling cream DAB	281140	1
Caelo	Cooling cream DAB	281822	2
Bombastus	Cooling cream DAB	281823	1
ichthyol/AHD	Cooling cream DAB	281921	2
Bomba/AHD	Cooling cream DAB	281920	1
Bomba/AHD	Cooling cream DAB	283115	1
Bombastus	Cooling cream DAB	284265	1
Anzag	Cooling cream DAB	5253700	1
Henry Lamotte/N/19.11.12	Cooling cream DAB	5253700	1
Anzag	Cooling cream DAB	60082802	1
Gehe	Cooling cream DAB 6 (stabili...	14268301	1
Gehe	Cooling cream DAB 6 (stabili...	14184501	1

- 87927 spectra from 500 *Apo-Ident* customers from a total of 26 040 batches from a further 332 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cooling cream / Cooling cream with rose oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cooling cream / Cooling cream with rose oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	700	0	34 844
Type B	0	3042	30	90 802
Type C	0	3361	64	87 927

The substance/substance group *Cooling cream* / *Cooling cream with rose oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.1429 %)
Type B	100.0000 % (> 99.9761 %)	99.0234 % (> 98.9258 %)
Type C	100.0000 % (> 99.6586 %)	98.1314 % (> 98.0438 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30889	33779	8.06	36.71
30890	31249	9.76	49.07
30921	31249	5.18	42.74
31052	31922	6.98	39.87
31224	31992	7.48	40.11
31249	31249	0.00	44.17
31457	33779	11.39	37.15
31471	31471	0.00	37.27
31922	31922	0.00	43.08
31992	31992	0.00	41.18
32294	32294	0.00	36.35
32506	32506	0.00	45.76
33779	33779	0.00	37.45

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Coriander oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31192-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Coriander oil; Oleum coriandri

Special notes

When selecting the *Coriander oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Coriander oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Coriander oil	12228509	31192	40	1407274
Caelo	Coriander oil	13288109	31407	40	20130903
Caelo	Coriander oil	170504	33913	40	20170228*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Coriander oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 161 spectra from 4 reference samples of the substance/substance group *Coriander oil*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Coriander oil	09127012	30322	41
Caelo	Coriander oil	10293303	30329	40
Caelo	Coriander oil	16140410	33253	40
Caelo	Coriander oil	181257	34056	40

- 93 713 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 77 spectra from 27 *Apo-Ident* customers from 40 batches from the substance/substance group *Coriander oil*.
- These include spectra of independent samples from 38 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Coriander oil	12228503	2
Caelo	Coriander oil	12061413	1
Caelo	Coriander oil	12228509	1
Caelo	Coriander oil	13288109	4
Caelo	Coriander oil	13288112	4
Caelo	Coriander oil	13288113	1
Caelo	Coriander oil	14319902	3
Caelo	Coriander oil	14319904	4
Caelo	Coriander oil	15229601	5
Caelo	Coriander oil	15229607	3
Caelo	Coriander oil	152229601	1
Taoasis	Coriander oil	104709-125892	1
Caelo	Coriander oil	15229608	1
Caelo	Coriander oil	1612210C/16140401	1
Caelo	Coriander oil	15229613	1
Caelo	Coriander oil	16140401	2
Caelo	Coriander oil	16140403	1
Caelo	Coriander oil	16140405	6
GEHE	Coriander oil	16140405	1
Caelo	Coriander oil	16140406	1
Caelo	Coriander oil	16140410	2
Caelo	Coriander oil	17050408	4
Caelo	Coriander oil	1711200C/17050403	1
Caelo	Coriander oil	2081306	1
Caelo	Coriander oil	7021303	1
Caelo	Coriander oil	10293307	2
Caelo	Coriander oil	11324003	1
Caelo	Coriander oil	17364701	1
Caelo	Coriander oil	11324002	1
Caelo	Coriander oil	11324004	1
Noweda	Coriander oil	11324005	1
Caelo	Coriander oil	11324005	4
AHD, 19.07.2013	Coriander oil	12228507	1
Caelo	Coriander oil	12228501	1
AHD, 26.09.2013	Coriander oil	12228509	2
AHD, 05.11.2013	Coriander oil	12228510	1
Kaders	Coriander oil	131528	1
Caelo	Coriander oil	14011403	1
Caelo	Coriander oil	18061201	1
Taoasis	Coriander oil	32040-108625	1
Taoasis	Coriander oil	31920-1111173	1
Caelo	Coriander oil	31081203	1
Taoasis	Coriander oil	33485-112618	1
VDI	Coriander oil	14319902	1

- 91 275 spectra from 501 *Apo-Ident* customers from a total of 26 969 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Coriander oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Coriander oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	161	0	93 713
Type C	0	75	2	91 275

The substance/substance group *Coriander oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.2733 %)
Type C	100.0000 % (> 99.6587 %)	97.4026 % (> 93.5065 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31192	31192	0.00	72.57
31407	33913	5.03	65.26
33913	33913	0.00	65.45

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **D-alpha tocopherol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30485-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

D-alpha tocopherol; alpha-Tocopherolum

Special notes

When selecting the *D-alpha tocopherol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *D-alpha tocopherol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	D-alpha tocopherol	13248808	31501	39	1507770
Caelo	D-alpha tocopherol	16005002	32072	40	20160203
Fagron	D-alpha tocopherol	16D14-B05-324514	32725	40	20160422
Fagron	D-alpha tocopherol	16H29-B08-334064	33419	40	20161027

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 159 spectra of 4 reference samples from the substance/substance group D-*alpha* tocopherol. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 385 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 160 spectra from 5 reference samples of the substance/substance group D-*alpha* tocopherol.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	D-alpha tocopherol	11313309	30485	40
Caelo	D-alpha tocopherol	13248808	31501 [†]	20
Caelo	D-alpha tocopherol	16005002	32072 [†]	20
Fagron	D-alpha tocopherol	16H29-B08-334064	33418	40
Fagron	D-alpha tocopherol	17B14-B02-336290	34036	40

- 93 714 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 60 spectra from 38 *Apo-Ident* customers from 37 batches from the substance/substance group *D-alpha tocopherol*.
- These include spectra of independent samples from 34 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	D-alpha tocopherol	13248804	1
Caelo	D-alpha tocopherol	16005002	9
Caelo	D-alpha tocopherol	16187604	4
Caelo	D-alpha tocopherol	8705007	1
Caelo	D-alpha tocopherol	11313308	3
Caelo	D-alpha tocopherol	712015	1
Sanacorp	D-alpha tocopherol	13248801	1
Caelo	D-alpha tocopherol	11313311	2
Caelo	D-alpha tocopherol	14159302	1
Caelo	D-alpha tocopherol	13248803	1
Caelo	D-alpha tocopherol	14159305	1
Caelo	D-alpha tocopherol	15102611	1
Caelo	D-alpha tocopherol	15102606	2
Gehe	D-alpha tocopherol	16187606	1
Audor Pharma	D-alpha tocopherol	16D14-B05-322951	1
Fagron	D-alpha tocopherol	16D14-B05-324514	1
Caelo	D-alpha tocopherol	17057501	1
Caelo	D-alpha tocopherol	17057508	1
Fagron	D-alpha tocopherol	17J31-B01	1
Caelo	D-alpha tocopherol	12106907	2
Frago/Noweda	D-alpha tocopherol	17B14-B02	1
Caelo	D-alpha tocopherol	12106901	1
Fagron	D-alpha tocopherol	20170816-03257	1
Caelo	D-alpha tocopherol	12106908	1
Caelo	D-alpha tocopherol	12106910	3
Caelo	D-alpha tocopherol	13248802	2
Caelo	D-alpha tocopherol	11313302	2
Caelo	D-alpha tocopherol	132488003	1
Caelo	D-alpha tocopherol	11313301	1
Caelo	D-alpha tocopherol	14159303	1
Caelo	D-alpha tocopherol	5012I-03257	1
Caelo	D-alpha tocopherol	13248801	1
Caelo	D-alpha tocopherol	14353802	1
Caelo	D-alpha tocopherol	14353806	2
Caelo	D-alpha tocopherol	16187606	1
Alte Apotheke	D-alpha tocopherol	16D14B05324514	1
Caelo	D-alpha tocopherol	18068402	1
Caelo	D-alpha tocopherol	11313309	1
Caelo	D-alpha tocopherol		1

- 91 292 spectra from 501 *Apo-Ident* customers from a total of 26 973 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are docu-

mented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *D-alpha tocopherol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *D-alpha tocopherol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	159	0	35 385
Type B	0	160	0	93 714
Type C	0	60	0	91 292

The substance/substance group *D-alpha tocopherol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2264 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.2500 %)
Type C	100.0000 % (> 99.6587 %)	100.0000 % (> 90.0000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31501	31501	0.00	55.05
32072	31501	6.86	56.68
32725	31501	17.84	70.32
33419	31501	29.08	76.66

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Dexpanthenol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30884-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Dexpanthenol; Dexpanthenolum

Special notes

When selecting the *Dexpanthenol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dexpanthenol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Dexpanthenol	12018707	30884	40	20120207
Caelo	Dexpanthenol	14051903	31317	40	20140319
Euro OTC	Dexpanthenol	1309030-03	31532	40	AR-15-FG-005869-01
Euro OTC	Dexpanthenol	1510008	31769	40	20151102*
Fagron	Dexpanthenol	17B06-B04-332078	33239	40	20170315

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Dexpanthenol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 344 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 477 spectra from 13 reference samples of the substance/substance group *Dexpanthenol*.
- These include spectra of independent samples from 11 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Dexpanthenol	1309030-03	31532 [†]	20
Euro OTC	Dexpanthenol	1510008	31769 [†]	17
Euro OTC	Dexpanthenol	1606015	32465	40
Euro OTC	Dexpanthenol	1607027	32467	40
Euro OTC	Dexpanthenol	1610041	32829	40
Euro OTC	Dexpanthenol	1703029	33238	40
Euro OTC	Dexpanthenol	1708013	33664	40
Caelo	Dexpanthenol	171721	33716	40
Caelo	Dexpanthenol	172877	33903	40
Euro OTC	Dexpanthenol	1801011	33929	40
Euro OTC	Dexpanthenol	1803008	34090	40
Euro OTC	Dexpanthenol	1806019-01	34210	40
Euro OTC	Dexpanthenol	1808007	34217	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 93 397 spectra from a total of 2154 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1383 spectra from 326 *Apo-Ident* customers from 384 batches from the substance/substance group *Dexpanthenol*.
- These include spectra of independent samples from 375 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Dexpanthenol	1304023	1
Caelo	Dexpanthenol	1302009	1
Euro OTC	Dexpanthenol	2081403	1
Caelo	Dexpanthenol	3815464	1
Caelo	Dexpanthenol	4031502	1
Euro OTC	Dexpanthenol	5051504	1
Euro OTC	Dexpanthenol	5051505	1
Euro OTC	Dexpanthenol	4091514	1
Euro OTC	Dexpanthenol	7011503	1
Euro OTC	Dexpanthenol	10061508	1
Euro OTC	Dexpanthenol	11021503	1
Euro OTC	Dexpanthenol	1111406	1
Caelo	Dexpanthenol	12031602	1
Caelo	Dexpanthenol	1235899	1
Caelo	Dexpanthenol	12285214	6
Alliance health	Dexpanthenol	1204010-01	1
Fagron	Dexpanthenol	12E11-N02	22
Fagron	Dexpanthenol	13/04-N04	3
Fagron	Dexpanthenol	13/29-B05-287528	1
Fagron	Dexpanthenol	12J30-B03-287502	1
Euro OTC	Dexpanthenol	1301009-01	24
Euro OTC	Dexpanthenol	130100901	2
Euro OTC	Dexpanthenol	1301009-02	12
Euro OTC	Dexpanthenol	130100902	2
Caelo	Dexpanthenol	13022-N12	1
Caelo	Dexpanthenol	13084709	6
Caelo	Dexpanthenol	13084703	7
Caelo	Dexpanthenol	13084705	3
Euro OTC	Dexpanthenol	130903001	2
Caesar/Noweda	Dexpanthenol	13084709	2
Phönix	Dexpanthenol	1309030-01	1
Euro OTC	Dexpanthenol	1309030-01	7
BASF AG	Dexpanthenol	1309030-01	1
Fagron	Dexpanthenol	1309030-02	1
EuRho	Dexpanthenol	1309030-01	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Dexpanthenol	1309030-02	10
Euro OTC	Dexpanthenol	130903002	1
Raben-Apotheke Mildenau	Dexpanthenol	1309030-03	1
Euro OTC	Dexpanthenol	1309030-03	3
Alliance Health 03.12.2...	Dexpanthenol	1309030-03	1
Caelo	Dexpanthenol	13189506	8
Caelo	Dexpanthenol	13189503	12
Caelo	Dexpanthenol	13342601	5
Caelo	Dexpanthenol	13342602	4
Caelo	Dexpanthenol	13342603	4
Caelo	Dexpanthenol	13342605	8
Caelo	Dexpanthenol	13342604	1
Bombastus	Dexpanthenol	13342605	1
Noweda	Dexpanthenol	13342605	1
Caelo	Dexpanthenol	13342606	12
Noweda	Dexpanthenol	13342606	1
Caelo	Dexpanthenol	13342610	6
Caelo	Dexpanthenol	13342608	1
GEHE 17.06.2014	Dexpanthenol	13342610	1
Caelo	Dexpanthenol	13342611	4
Fagron	Dexpanthenol	13D11-N12	1
Fagron	Dexpanthenol	13D22-N12	14
Fagron	Dexpanthenol	13I04-N04	15
Fagron	Dexpanthenol	13i04-n04	1
Fagron	Dexpanthenol	13J04-N04	2
Fagron	Dexpanthenol	13J29-B05	8
Fagron	Dexpanthenol	13j29-B05-287525	1
Fagron	Dexpanthenol	13J29-B05-287526	7
Caelo	Dexpanthenol	13J29-B05-287526	1
Fagron	Dexpanthenol	13J29-b05-287526	1
Fagron	Dexpanthenol	13J29B05287526	1
Anzag	Dexpanthenol	13j29-B05-287527	1
	Dexpanthenol	13J29-B05-287527	6
Fagron	Dexpanthenol	13J29-B05-287527	2
	Dexpanthenol	13J29B05287527	1
Fagron	Dexpanthenol	13J29-B05-287528	1
Fagron	Dexpanthenol	13L02-B04	3
Fagron	Dexpanthenol	13L02-B04-292067	2
Fagron	Dexpanthenol	13L02-B04-292068	8
Fagron	Dexpanthenol	13L02-BB04-292067	1
Caelo	Dexpanthenol	14031410	1
Caelo	Dexpanthenol	14051901	1
Caelo	Dexpanthenol	14051903	6
Caelo	Dexpanthenol	14051904	6
Caelo	Dexpanthenol	14051905	7
Fagron	Dexpanthenol	14051908	2
Caelo	Dexpanthenol	14051908	1
Caelo	Dexpanthenol	14051913	1
Caelo	Dexpanthenol	14051909	2
Euro OTC	Dexpanthenol	14070390	1
Euro OTC	Dexpanthenol	1407039-01	10
Euro OTC	Dexpanthenol	140703901	3
Caelo	Dexpanthenol	1407039-01	2
BASF	Dexpanthenol	1407039-01	1
Euro OTC	Dexpanthenol	140703902	1
Nachprüfung	Dexpanthenol	1407039-01	1
	Dexpanthenol	1407039-02	2
Euro OTC	Dexpanthenol	1407039-02	7
Alliance Health 27.03.2...	Dexpanthenol	1407039-02	1
Sanacorp	Dexpanthenol	1407039-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
EuRho/Noweda	Dexpanthenol	1407039-02	1
Euro OTC	Dexpanthenol	1407039-1	1
Euro OTC	Dexpanthenol	140739	1
Euro OTC	Dexpanthenol	14070399-01	1
Caelo	Dexpanthenol	14214501	3
Caelo	Dexpanthenol	14214502	5
Caelo	Dexpanthenol	142141508	1
Caelo	Dexpanthenol	14214504	6
Caelo	Dexpanthenol	14214505	6
Caelo	Dexpanthenol	14214507	3
Caelo	Dexpanthenol	14214508	2
Caelo	Dexpanthenol	14214510	3
Fagron	Dexpanthenol	14214510	1
Caelo	Dexpanthenol	14214511	5
Caelo	Dexpanthenol	14214512	4
Phoenix goettin	Dexpanthenol	14214512	1
Caelo	Dexpanthenol	1421508	1
14356201	Dexpanthenol	14356201	2
Caelo	Dexpanthenol	14356201	8
Caelo	Dexpanthenol	14356203	2
Caelo	Dexpanthenol	14356204	4
Caelo	Dexpanthenol	14356210	3
Fagron	Dexpanthenol	14E01-811-295577	1
Fagron	Dexpanthenol	14e01-811-300068	1
Fagron	Dexpanthenol	14E01-B11	8
Caelo	Dexpanthenol	14E01-B11	1
Sanacorp	Dexpanthenol	14E01-B11	1
Fagron	Dexpanthenol	14e01-b11-293623	1
anzag	Dexpanthenol	14E01-B11-293624	7
Fagron	Dexpanthenol	14E01B11293623	1
Fagron	Dexpanthenol	14E01-B11-293623	1
Fagron	Dexpanthenol	14E01-B11-293624	2
Fagron	Dexpanthenol	14E01-B11-295534	1
Fagron	Dexpanthenol	14E01-B11-295533	1
Anzag	Dexpanthenol	14E01-B11-295534	3
Fagron	Dexpanthenol	14E01-B11-295544	2
Fagron	Dexpanthenol	14E01-B11-295577	2
Fagron	Dexpanthenol	14E01-b11-295544	1
Fagron	Dexpanthenol	14E01-B11-295578	10
Fargon/Alliance HC	Dexpanthenol	14E01-B11295578	1
Fagron	Dexpanthenol	14E01-B11-3000068	1
Fagron	Dexpanthenol	14E01-B11-300065	2
A.H.D.	Dexpanthenol	14E01-B11-300068	1
Fagron	Dexpanthenol	14E01-B11-300068	12
Caelo	Dexpanthenol	14H05-B02-3000066	1
Fagron	Dexpanthenol	14H05-B02-0001	1
apomix/Holdermann	Dexpanthenol	14H05-B02-300066	1
Fagron	Dexpanthenol	14H05-B02	2
Fagron	Dexpanthenol	14Eo1-B11-300068	1
Caelo	Dexpanthenol	14H05B02300069	1
Fagron	Dexpanthenol	14H05-B02-300066	1
Anzag	Dexpanthenol	14H05-B02-300067	1
Fagron	Dexpanthenol	14H05-B02-300069	11
Caelo	Dexpanthenol	14H05-B02-300069	1
Anzag	Dexpanthenol	14H05-B02-305656	3
Fagron	Dexpanthenol	14H05-B02-305656	2
Caelo	Dexpanthenol	14H05-B02-305656	1
Euro OTC	Dexpanthenol	1502010-01	1
Euro OTC	Dexpanthenol	1502010-02	1
Euro OTC	Dexpanthenol	1502019-01	20

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe/Euro 30.06.15	Dexpanthenol	1502019-01	1
Caelo	Dexpanthenol	1502019-01	2
Phönix Gotha	Dexpanthenol	1502019-01	1
Euro OTC	Dexpanthenol	150201901	1
Euro OTC	Dexpanthenol	1502019-02	17
Phönix Gotha	Dexpanthenol	1502019-02	1
Euro OTC	Dexpanthenol	150201902	2
BASF/ Phönix	Dexpanthenol	1502019-02	2
Bombastus	Dexpanthenol	1502019-02	2
Euro OTC	Dexpanthenol	150219-02	1
Caelo	Dexpanthenol	15097701	6
Caelo	Dexpanthenol	15097703	7
Caelo	Dexpanthenol	15097704	7
Caelo	Dexpanthenol	15097705	5
Caelo	Dexpanthenol	15097708	8
Fagron	Dexpanthenol	15097709	1
Caelo	Dexpanthenol	15097709	3
Caelo	Dexpanthenol	15097710	5
Caelo	Dexpanthenol	15097711	9
Caelo	Dexpanthenol	15097712	2
Noweda	Dexpanthenol	15097712	1
Caelo	Dexpanthenol	15097713	3
Caelo	Dexpanthenol	15097714	13
Caesar & Loretz GmbH Hi...	Dexpanthenol	15097714	2
Sanacorp	Dexpanthenol	15097714	1
Caelo	Dexpanthenol	15097716	1
Euro OTC	Dexpanthenol	15100008-01	1
Euro OTC	Dexpanthenol	15100008-02	2
Noweda	Dexpanthenol	1510008	1
Euro OTC	Dexpanthenol	1510008-01	27
Euro OTC	Dexpanthenol	151000801	1
Caelo	Dexpanthenol	1510008-01	2
Gehe/ Euro	Dexpanthenol	1510008-01	1
BASF AG/Noweda	Dexpanthenol	1510008-01	1
Euro OTC	Dexpanthenol	1510008-02	13
Euro OTC	Dexpanthenol	151000802	1
Phoenix	Dexpanthenol	1510008-02	1
Euro OTC	Dexpanthenol	151008-01	2
Caelo	Dexpanthenol	15260905	6
Alliance H.	Dexpanthenol	15260905	1
Caelo	Dexpanthenol	15260906	10
Caelo	Dexpanthenol	15260907	5
Caelo	Dexpanthenol	15260908	1
Euro OTC	Dexpanthenol	15260908	1
Caelo	Dexpanthenol	15395502	10
Caelo	Dexpanthenol	15395504	3
Caolo	Dexpanthenol	15395506	1
Caelo	Dexpanthenol	15395506	1
Caelo	Dexpanthenol	15395507	3
Caelo	Dexpanthenol	15395508	6
Fagron	Dexpanthenol	15A30-B08	4
phönix	Dexpanthenol	15a30b08312353	1
Fagron	Dexpanthenol	15A30-B08-306694	4
Fagron	Dexpanthenol	15A30-B08-306695	6
Fagron	Dexpanthenol	15A30-B08-310962	2
Alliance	Dexpanthenol	15A30-B08-310962	2
Fagron	Dexpanthenol	15A30-B08-312353	5
Caelo	Dexpanthenol	15A30-B08-312353	1
Fagron	Dexpanthenol	15F26-B01-314370	9
Euro OTC	Dexpanthenol	15F26-B01-314370	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
DSM Nutritional Products UK	Dexpanthenol	15F26-B01-314370	1
Fagron	Dexpanthenol	15F26-B01-317519	2
Alliance	Dexpanthenol	15F26-B01-317520	2
Fagron	Dexpanthenol	15F26-B01-317520	1
Fagron	Dexpanthenol	15F26-B01-319991	6
Caelo	Dexpanthenol	1604305	1
Caelo	Dexpanthenol	16043502	5
Caelo	Dexpanthenol	16043503	3
Caelo	Dexpanthenol	16043504	3
Caelo	Dexpanthenol	16043505	4
Noweda	Dexpanthenol	16043505	1
DSM Nutritional Prod./...	Dexpanthenol	16043505	1
Hedinger	Dexpanthenol	16043505	1
Caelo	Dexpanthenol	16043510	9
Caelo	Dexpanthenol	16043506	1
Caelo	Dexpanthenol	16043511	3
Caelo	Dexpanthenol	16043512	9
caesar	Dexpanthenol	16043512	1
Euro OTC	Dexpanthenol	1606015-01	17
BASF AG/Noweda	Dexpanthenol	1607027-01	1
Euro OTC	Dexpanthenol	1607027-01	11
Euro OTC	Dexpanthenol	160702701	2
Euro/Noweda30.01.17	Dexpanthenol	1606015-01	1
GEHE	Dexpanthenol	1607027-01	1
Caelo	Dexpanthenol	1610041-01	4
Euro OTC	Dexpanthenol	161004101	2
Euro OTC	Dexpanthenol	1610041-01	14
hochhaus-Apotheke	Dexpanthenol	1610041-01	1
GEHE	Dexpanthenol	161004101	1
Euro OTC	Dexpanthenol	1610041-02	8
Euro OTC	Dexpanthenol	161220D/1606015-01	1
Caelo	Dexpanthenol	16164808	3
Caelo	Dexpanthenol	16164809	2
Caelo	Dexpanthenol	16164810	3
Caelo	Dexpanthenol	16164811	2
Caelo	Dexpanthenol	16164813	5
Caelo	Dexpanthenol	16164816	1
Caelo	Dexpanthenol	16311603	2
Caelo	Dexpanthenol	16311604	4
Caelo	Dexpanthenol	16311605	8
Caelo	Dexpanthenol	16311606	1
Caelo	Dexpanthenol	16311610	2
DSM Nutritional Prod./...	Dexpanthenol	16311610	1
Caelo	Dexpanthenol	16D07-804-322508	1
Fagron	Dexpanthenol	16D07-B04	1
Fagron	Dexpanthenol	16D07-B04-322110	4
Caelo	Dexpanthenol	16D07-B04-322110	1
phönix	Dexpanthenol	16D07-B04-322110	1
Fagron	Dexpanthenol	16D07-B04-322508	9
Noweda	Dexpanthenol	16D07B04322508	1
Fagron	Dexpanthenol	16D07-B04-323096	6
Alliance	Dexpanthenol	16D07-B04-323097	1
Fagron	Dexpanthenol	16D07-B04-323097	2
Fagron	Dexpanthenol	16D07-B04-327322	3
Sanacorp	Dexpanthenol	1703029-01	1
Caelo	Dexpanthenol	16Do07-B04-323096	1
Fagron	Dexpanthenol	17/21-B01-340266	1
Euro OTC	Dexpanthenol	1703029-01	7
Fagron	Dexpanthenol	17/21-B01-340264	1
Fagron	Dexpanthenol		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Dexpanthenol	1703029-01	1
beyvers	Dexpanthenol	1703029-01	1
apo ident	Dexpanthenol	1703029-01	1
Euro OTC	Dexpanthenol	17041503	1
Caelo	Dexpanthenol	17065403	2
Caelo	Dexpanthenol	17065404	2
Caelo	Dexpanthenol	17065406	3
Caelo	Dexpanthenol	17065405	1
Caelo	Dexpanthenol	17065411	2
Caelo	Dexpanthenol	17065412	6
DSM Nutritional Prod./...	Dexpanthenol	17065412	1
Euro OTC	Dexpanthenol	1708013-01	4
GEHE	Dexpanthenol	17065413	1
Caelo	Dexpanthenol	1708013-01	1
Euro OTC	Dexpanthenol	171106D/1703029-01	1
Fagron	Dexpanthenol	171120D/TL01601003	1
Euro OTC	Dexpanthenol	17091401	1
Caelo	Dexpanthenol	17172102	1
Caelo	Dexpanthenol	17172103	3
Caelo	Dexpanthenol	17172106	1
Caelo	Dexpanthenol	17172107	2
Caelo	Dexpanthenol	17172110	1
Phönix	Dexpanthenol	1721b01-340264	1
Caelo	Dexpanthenol	172587709	1
Caelo	Dexpanthenol	17287707	1
Caelo	Dexpanthenol	17335101	1
Caelo	Dexpanthenol	17335102	1
Caelo	Dexpanthenol	17335103	1
Gehe	Dexpanthenol	17335105	1
Gehe	Dexpanthenol	17335103	1
Fagron	Dexpanthenol	17B06-B04	3
Fagron	Dexpanthenol	17B0619	1
Caelo	Dexpanthenol	17335105	2
Noweda	Dexpanthenol	17B06-B04-322078	1
Fagron	Dexpanthenol	17B06-B04-332075	5
Sanaorp	Dexpanthenol	17B06-B04-332077	2
Fagron	Dexpanthenol	17B06-B04-332077	3
Caelo	Dexpanthenol	17B06-B04-332076	1
Fagron	Dexpanthenol	17B06-B04-332076	1
Fagron	Dexpanthenol	17B06-B04-332078	3
Noweda	Dexpanthenol	17B06-B04-332078	2
Fagron	Dexpanthenol	17I21-B01-340265	1
Noweda	Dexpanthenol	17D07-B04-322508	1
Noweda	Dexpanthenol	17I21B01340269	1
Fagron	Dexpanthenol	17I21-B01-340269	1
Fagron	Dexpanthenol	17I21-B01-340264	1
Euro OTC	Dexpanthenol	1801011-01	2
Sanacorp	Dexpanthenol	2/AKD4	1
BASF	Dexpanthenol	2181	1
Caelo	Dexpanthenol	2187	1
Caelo	Dexpanthenol	22011511	1
Caelo	Dexpanthenol	23041401	1
Euro OTC	Dexpanthenol	23041402	1
Caelo	Dexpanthenol	23081603	1
Caelo	Dexpanthenol	23081602	1
Euro OTC	Dexpanthenol	23101409	1
Caelo	Dexpanthenol	26081401	1
Caelo	Dexpanthenol	26081402	1
Caelo	Dexpanthenol	2614E-01082	1
Euro OTC	Dexpanthenol	27031503	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Dexpanthenol	28071501	1
Euro OTC	Dexpanthenol	28111404	1
Fagron	Dexpanthenol	2914E-01082	1
Fagron	Dexpanthenol	295544	1
Euro OTC	Dexpanthenol	31071501	1
BASF/Noweda	Dexpanthenol	5116I-01082	1
BASF AG	Dexpanthenol	94515888Q0	1
Caelo	Dexpanthenol		3
Sanacorp	Dexpanthenol		1
ACM	Dexpanthenol	KNC0684J15	1
ACM	Dexpanthenol	KNC1152J15	1
Fagron	Dexpanthenol	TL01203018	4
DSM Nutritional Prod./...	Dexpanthenol	TL01212052	1
Fagron	Dexpanthenol	TL01305030	8
Caelo	Dexpanthenol	TL01503022	1
Caelo	Dexpanthenol	TL01508051	1
Fagron	Dexpanthenol	TL01312070	1
Fagron	Dexpanthenol	TL01503022	1
Fagron	Dexpanthenol	TL01601003	1
Caelo	Dexpanthenol	3101202	1
sanaEuro	Dexpanthenol	311V-01082	1
Caelo	Dexpanthenol	3051305	1
Caelo	Dexpanthenol	2101303	1
Fagron	Dexpanthenol	TL01611112	1
Caelo	Dexpanthenol	3121303	1
Caelo	Dexpanthenol	5061307	1
Caelo	Dexpanthenol	5061306	1
Fagron	Dexpanthenol	8572983	1
Caelo	Dexpanthenol	9081206	1
Caelo	Dexpanthenol	9312203	1
Caelo	Dexpanthenol	10067908	1
Noweda	Dexpanthenol	11010703	1
Kehr	Dexpanthenol	110032-01	1
EuRho	Dexpanthenol	1110032-01	1
EuRho	Dexpanthenol	1107023-01	1
Sanacorp/ euRho	Dexpanthenol	1103053-01	1
Euro OTC	Dexpanthenol	1110032-01	3
Caelo	Dexpanthenol	11293501	2
Gehe	Dexpanthenol	11293501	1
Caelo	Dexpanthenol	11293506	3
Caelo	Dexpanthenol	11293503	1
Caelo	Dexpanthenol	11293507	4
Caelo	Dexpanthenol	1193507	1
Fagron	Dexpanthenol	11H26-N03	4
Fagron	Dexpanthenol	11B21-N01	1
Caelo	Dexpanthenol	12018702	1
Fagron	Dexpanthenol	11B21-N01	1
Caelo	Dexpanthenol	12018703	10
Caelo	Dexpanthenol	12018706	3
Caelo	Dexpanthenol	12018707	4
Caelo	Dexpanthenol	12018713	5
Caelo	Dexpanthenol	12018715	1
Caelo	Dexpanthenol	1201873	1
Euro OTC	Dexpanthenol	1204010-01	9
Phönix Gotha	Dexpanthenol	1204010-01	1
Euro OTC	Dexpanthenol	1204010-02	12
EuRho	Dexpanthenol	1204010-02	1
Caelo	Dexpanthenol	12126702	14
Caelo	Dexpanthenol	12126704	5
Caelo	Dexpanthenol	12126705	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Anzag	Dexpanthenol	12126705	2
Caelo	Dexpanthenol	12126707	2
Caelo	Dexpanthenol	1213A-01082	1
Großhandel	Dexpanthenol	12186705	1
Phönix	Dexpanthenol	1216707	1
Caelo	Dexpanthenol	12285201	6
Caesar&loretz,Phoenix	Dexpanthenol	12285201	1
Caelo	Dexpanthenol	12285202	7
Caelo	Dexpanthenol	12285203	4
Gehe	Dexpanthenol	12285203	1
cesar&lorentz/Phönix	Dexpanthenol	12285203	1
Caelo	Dexpanthenol	12285204	3
Phönix	Dexpanthenol	12285204	1
Caelo	Dexpanthenol	12285205	1
Noweda	Dexpanthenol	12285211	2
Caelo	Dexpanthenol	12285211	8
Caelo	Dexpanthenol	12285212	2
Caelo	Dexpanthenol	122852211	1
Euro OTC	Dexpanthenol	129401002	1
Fagron	Dexpanthenol	12813-No7	1
Fagron	Dexpanthenol	12B13-N07	4
Noweda	Dexpanthenol	12B13-N07	1
Fagron	Dexpanthenol	12b13-n07	1
Gehe	Dexpanthenol	12E11-N02	2
Fagron	Dexpanthenol	12e11-n02	1
Fragron	Dexpanthenol	12E11-N02	1
Anzag	Dexpanthenol	12E11-N02	3
fargon/Sanacorp	Dexpanthenol	12E11-N02	1
Noweda	Dexpanthenol	12E11-N02	3
Caelo	Dexpanthenol	12E11-N02	1
Fagron	Dexpanthenol	12e11-N02	1
Fagron	Dexpanthenol	12E11-No2	1
Fagron	Dexpanthenol	12E11-N02	1
Fagron	Dexpanthenol	12H28-N07	6
Phoenix 30.04.2013	Dexpanthenol	12H28-N07	1
Fagron	Dexpanthenol	12I04-N04	1
Anzag	Dexpanthenol	12J30-B03-287502	1
Fagron	Dexpanthenol	12K27-N04	7
Anzag	Dexpanthenol	12K27-N04	1
Fagron	Dexpanthenol	12L28-N02	1
Fagron	Dexpanthenol	12L28-N03	6
Caelo	Dexpanthenol	12L28-N03	1
Gehe	Dexpanthenol	12L28-N03	1
Fagron	Dexpanthenol	12L28-No3	1
Euro OTC	Dexpanthenol	1301009	1
Fagron	Dexpanthenol	12L28-N03	1
Anzag	Dexpanthenol	13/04-N04	1
Gehe	Dexpanthenol	13/04-N04	1
Euro OTC	Dexpanthenol	1301009-001	1
Alliance	Dexpanthenol	1301009-01	1
Daiichi Fine Chemicals	Dexpanthenol	1301009-01	1
EuRho	Dexpanthenol	1301009-01	1
Caelo	Dexpanthenol	1301009-02	1
Phönix	Dexpanthenol	1301009-02	1
Caelo	Dexpanthenol	13021305	1
Caelo	Dexpanthenol	13084702	15
Großhandel	Dexpanthenol	13084702	1
Caelo	Dexpanthenol	13084704	6
Europ Ots Pharma	Dexpanthenol	13084705	1
Caelo	Dexpanthenol	13084706	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dexpanthenol	13084707	2
Caelo	Dexpanthenol	13084711	1
Caelo	Dexpanthenol	1309030-03	1
Caelo	Dexpanthenol	13189502	2
Caelo	Dexpanthenol	13091309	1
Großhandel	Dexpanthenol	13189502	1
Großhandel	Dexpanthenol	13189503	2
Caelo	Dexpanthenol	13189504	2
Noweda	Dexpanthenol	13342601	1
Sanacorp	Dexpanthenol	13342601	1
Gehe	Dexpanthenol	13D22-N12	2
Caelo	Dexpanthenol	1513F-01082	1
Caelo	Dexpanthenol	15071304	1
Caelo	Dexpanthenol	16021307	1
Caelo	Dexpanthenol	16101305	1
Caelo	Dexpanthenol	18111306	1
Caelo	Dexpanthenol	18011303	1
Caelo	Dexpanthenol	19121204	1
Caelo	Dexpanthenol	20071205	1
Caelo	Dexpanthenol	20061203	1
Caelo	Dexpanthenol	19121205	1
Caelo	Dexpanthenol	20121301	1
Caelo	Dexpanthenol	20261012-1	1
Caelo	Dexpanthenol	22021304	1
Fagron	Dexpanthenol	246019-0001	1
Caelo	Dexpanthenol	25021405	1
Caelo	Dexpanthenol	24111201	1
Caelo	Dexpanthenol	25101204	1
Caelo	Dexpanthenol	25031305	2
Sanacorp	Dexpanthenol	36/AKD22	1
Caelo	Dexpanthenol	25101203	1
Caelo	Dexpanthenol	3713A-01082	1
Euroo/Ebert	Dexpanthenol	3813E-01082	1
Caesar/Noweda	Dexpanthenol	4811A-01082	1
Caelo	Dexpanthenol	82729498	1
Euro OTC	Dexpanthenol	E0EIE22902	1
Euro OTC	Dexpanthenol	95700556P0	2
Euro OTC	Dexpanthenol	E0EJD41902	1
Gehe	Dexpanthenol	TL01108059	1
Fagron	Dexpanthenol	TL01110074	6
Caelo	Dexpanthenol	TL01110076	1
Fagron	Dexpanthenol	T1012120252	1
Fagron	Dexpanthenol	TL01110074	1
Fagron	Dexpanthenol	TL01212052	1
Fagron	Dexpanthenol	TL01212052	1
Anzag	Dexpanthenol	12H28-N07	1
Fagron	Dexpanthenol	12e11-no2	1
Fagron	Dexpanthenol	12K27-N04	1
Holdermann	Dexpanthenol	13/04-N04	1
EUro/AHD	Dexpanthenol	1301009-01	2
Hedinger/AHD	Dexpanthenol	1301009-02	1
A.H.D.	Dexpanthenol	13D22-N12	1
Phönix	Dexpanthenol	13J29-B05287525	1

- 89 969 spectra from 501 *Apo-Ident* customers from a total of 26 630 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dexpantenol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dexpantenol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	477	0	93 397
Type C	0	1371	12	89 969

The substance/substance group *Dexpantenol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.7421 %)
Type C	100.0000 % (> 99.6586 %)	99.1323 % (> 98.9154 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30884	31532	3.77	129.01
31317	31532	1.64	139.15
31532	31532	0.00	142.43
31769	31769	0.00	129.81
33239	31532	1.95	139.44

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Dexpanthenol concentrate 50%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30492-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Dexpanthenol concentrate 50%; Dexpanthenolum 50%; Pantothenolum 50%

Special notes

When selecting the *Dexpanthenol concentrate 50%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dexpanthenol concentrate 50%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Dexpanthenol con...	12171803	30664	40	20120724
Caelo	Dexpanthenol con...	15408002	31835	40	20151211
Caelo	Dexpanthenol con...	153314	31879	40	20151014*
Caelo	Dexpanthenol con...	162483	32579	40	20160901*
Caelo	Dexpanthenol con...	18050802	34130	40	20180213

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Dexpanthenol concentrate 50%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 344 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 324 spectra from 9 reference samples of the substance/substance group *Dexpanthenol concentrate 50%*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dexpanthenol concentrate 50%	11358202	30408	40
Caelo	Dexpanthenol concentrate 50%	11279706	30492	40
Caelo	Dexpanthenol concentrate 50%	14116302	31461	60
Caelo	Dexpanthenol concentrate 50%	15408002	31835 [†]	20
Caelo	Dexpanthenol concentrate 50%	153314	31879 [†]	20
Caelo	Dexpanthenol concentrate 50%	163336	33161	40
Caelo	Dexpanthenol concentrate 50%	173261	33974	40
Caelo	Dexpanthenol concentrate 50%	172781	33975	40
Caelo	Dexpanthenol concentrate 50%	16248301	not required	24

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 550 spectra from a total of 2158 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 474 spectra from 119 *Apo-Ident* customers from 85 batches from the substance/substance group *Dexpanthenol concentrate 50%*.
- These include spectra of independent samples from 80 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Dexpanthenol concentrate 50%	2101404	1
Caelo	Dexpanthenol concentrate 50%	11358201	2
Caelo	Dexpanthenol concentrate 50%	12081406	1
Nachprüfung	Dexpanthenol concentrate 50%	13019203	2
Caelo	Dexpanthenol concentrate 50%	13100408	1
Caelo	Dexpanthenol concentrate 50%	13305003	10
Caelo	Dexpanthenol concentrate 50%	13382102	5
Gehe 17.06.2014	Dexpanthenol concentrate 50%	13382102	1
Caelo	Dexpanthenol concentrate 50%	14020803	3
Anzag	Dexpanthenol concentrate 50%	14020803	2
Caelo	Dexpanthenol concentrate 50%	1402003	1
Euro OTC	Dexpanthenol concentrate 50%	14020803	2
VDL	Dexpanthenol concentrate 50%	14020804	1
Caelo	Dexpanthenol concentrate 50%	14020804	7
sanacorp	Dexpanthenol concentrate 50%	14116301	1
Caelo	Dexpanthenol concentrate 50%	14116301	6
Sanacorp	Dexpanthenol concentrate 50%	14116302	1
Caelo	Dexpanthenol concentrate 50%	14116302	6
Raben-Apotheke Mildenau	Dexpanthenol concentrate 50%	14116302	1
Caelo	Dexpanthenol concentrate 50%	14198101	14
Caelo	Dexpanthenol concentrate 50%	14198103	4
Fagron	Dexpanthenol concentrate 50%	14198103	1
sanacorp	Dexpanthenol concentrate 50%	14198103	1
Caelo	Dexpanthenol concentrate 50%	14199101	1
Caelo	Dexpanthenol concentrate 50%	14275203	9
Caelo	Dexpanthenol concentrate 50%	14275202	1
Caelo	Dexpanthenol concentrate 50%	14275204	4
Caelo	Dexpanthenol concentrate 50%	14275207	4
Caelo	Dexpanthenol concentrate 50%	14373301	5
Caelo	Dexpanthenol concentrate 50%	14373302	7
Gehe	Dexpanthenol concentrate 50%	14373302	1
Caelo	Dexpanthenol concentrate 50%	15050802	20
Caelo	Dexpanthenol concentrate 50%	15050801	2
Caelo	Dexpanthenol concentrate 50%	1508002	1
Caelo	Dexpanthenol concentrate 50%	15142901	6
Caelo	Dexpanthenol concentrate 50%	15142902	3
Euro OTC	Dexpanthenol concentrate 50%	15142902	1
Caelo	Dexpanthenol concentrate 50%	15196202	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dexpanthenol concentrate 50%	1514902	1
Caelo	Dexpanthenol concentrate 50%	15196203	14
Noweda	Dexpanthenol concentrate 50%	15196202	1
Caelo	Dexpanthenol concentrate 50%	15331402	2
Caelo	Dexpanthenol concentrate 50%	15331404	14
Rondell Apotheke	Dexpanthenol concentrate 50%	15331404	1
Caelo	Dexpanthenol concentrate 50%	15408002	13
Caesar&Loretz GmbH,Phön...	Dexpanthenol concentrate 50%	15408002	1
Noweda	Dexpanthenol concentrate 50%	15408002	1
Caelo	Dexpanthenol concentrate 50%	15408003	9
Caesar&Loretz GmbH	Dexpanthenol concentrate 50%	15408003	1
Caelo	Dexpanthenol concentrate 50%	1602503	1
Caelo	Dexpanthenol concentrate 50%	16025303	9
Bombastus	Dexpanthenol concentrate 50%	16025303	1
Ebert & Jacobi 01.06...	Dexpanthenol concentrate 50%	16025303	1
Caelo	Dexpanthenol concentrate 50%	16136503	16
Caelo	Dexpanthenol concentrate 50%	16136504	15
Caelo	Dexpanthenol concentrate 50%	16136505	7
Sanacorp	Dexpanthenol concentrate 50%	16136505	1
caesar&Loretz GmbH,AHC,...	Dexpanthenol concentrate 50%	16248301	1
Caesar & Loretz GmbH	Dexpanthenol concentrate 50%	16248301	1
Caelo	Dexpanthenol concentrate 50%	16248301	7
Caelo	Dexpanthenol concentrate 50%	16248303	18
Caesar & Loretz GmbH	Dexpanthenol concentrate 50%	16248303	1
Caesar & Loretz GmbH	Dexpanthenol concentrate 50%	16315102	2
Caelo	Dexpanthenol concentrate 50%	16315102	5
Caelo	Dexpanthenol concentrate 50%	16333602	9
Caelo	Dexpanthenol concentrate 50%	16333603	23
Fagron	Dexpanthenol concentrate 50%	16333603	1
Caesar & Loretz GmbH	Dexpanthenol concentrate 50%	16333603	1
Caelo	Dexpanthenol concentrate 50%	1633603	1
Caelo	Dexpanthenol concentrate 50%	17031301	4
Caelo	Dexpanthenol concentrate 50%	17059201	1
Caelo	Dexpanthenol concentrate 50%	17059202	10
Finteler Apotheke	Dexpanthenol concentrate 50%	17031303	1
Caelo	Dexpanthenol concentrate 50%	17031303	2
Caesar&Loretz GmbH,27.0...	Dexpanthenol concentrate 50%	17059203	1
Caelo	Dexpanthenol concentrate 50%	17059203	4
Caesar & Loretz GmbH	Dexpanthenol concentrate 50%	17135701	1
Caelo	Dexpanthenol concentrate 50%	17135701	1
Caesar & Loretz GmbH	Dexpanthenol concentrate 50%	17135702	1
Caelo	Dexpanthenol concentrate 50%	17135702	9
Caelo	Dexpanthenol concentrate 50%	17194601	6
Alliance	Dexpanthenol concentrate 50%	17194601	1
Caelo	Dexpanthenol concentrate 50%	17194602	2
Fagron	Dexpanthenol concentrate 50%	17194602	1
Caelo	Dexpanthenol concentrate 50%	17194603	3
Caelo	Dexpanthenol concentrate 50%	17278101	4
Caelo	Dexpanthenol concentrate 50%	17278103	1
Caelo	Dexpanthenol concentrate 50%	17326101	2
Caelo	Dexpanthenol concentrate 50%	17326102	5
GEHE	Dexpanthenol concentrate 50%	17326102	1
Caelo	Dexpanthenol concentrate 50%	18050802	1
Audor Pharma	Dexpanthenol concentrate 50%	AP65524488Q0P	2
Lamotte	Dexpanthenol concentrate 50%	AP83943088Q0P	1
Audor Pharma	Dexpanthenol concentrate 50%	AP36716516K0P	1
Sanacorp	Dexpanthenol concentrate 50%		1
Audor Pharma	Dexpanthenol concentrate 50%	AP83943088Q0P	1
Caelo	Dexpanthenol concentrate 50%		1
Caelo	Dexpanthenol concentrate 50%	11279703	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dexpanthenol concentrate 50%	11279704	1
Caelo	Dexpanthenol concentrate 50%	11358202	1
Gehe	Dexpanthenol concentrate 50%	11358204	1
Caelo	Dexpanthenol concentrate 50%	12030801	2
Caelo	Dexpanthenol concentrate 50%	12127003	4
Caelo	Dexpanthenol concentrate 50%	12171802	3
Noweda	Dexpanthenol concentrate 50%	12171802	1
Caelo	Dexpanthenol concentrate 50%	12171803	3
Caelo	Dexpanthenol concentrate 50%	12256001	4
Kehr	Dexpanthenol concentrate 50%	12171803	2
Caelo	Dexpanthenol concentrate 50%	12256002	3
Alliance Healthcare	Dexpanthenol concentrate 50%	12256001	1
Caelo	Dexpanthenol concentrate 50%	12340302	4
Noweda	Dexpanthenol concentrate 50%	12340302	1
Sanacorp	Dexpanthenol concentrate 50%	12340302	1
Noweda	Dexpanthenol concentrate 50%	12340303	2
Caelo	Dexpanthenol concentrate 50%	13019203	7
Gehe	Dexpanthenol concentrate 50%	13019203	1
Caelo	Dexpanthenol concentrate 50%	13019204	1
Gehe	Dexpanthenol concentrate 50%	13019204	1
Caelo	Dexpanthenol concentrate 50%	13068003	5
Alliance Health	Dexpanthenol concentrate 50%	13068003	1
Caelo	Dexpanthenol concentrate 50%	13068004	4
Caelo	Dexpanthenol concentrate 50%	13174702	7
gehe 31.01.2014	Dexpanthenol concentrate 50%	13068004	1
Caelo	Dexpanthenol concentrate 50%	13174703	7
Alliance Healthcare	Dexpanthenol concentrate 50%	13174703	1
Sonnenberg-Apotheke	Dexpanthenol concentrate 50%	13305004	1
Caelo	Dexpanthenol concentrate 50%	13305004	4
Caesar&loretz,Phoenix	Dexpanthenol concentrate 50%	14116302	2
Sanacorp	Dexpanthenol concentrate 50%	AKD29	1
Caelo	Dexpanthenol concentrate 50%	D-050413-2	1
Audor Pharma	Dexpanthenol concentrate 50%	AP83943088QOP	1

- 90 878 spectra from 501 *Apo-Ident* customers from a total of 26 925 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dexpanthenol concentrate 50%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dexpanthenol concentrate 50%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	324	0	93 550
Type C	0	473	1	90 878

The substance/substance group *Dexpanthenol concentrate 50%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.1481 %)
Type C	100.0000 % (> 99.6586 %)	99.7890 % (> 99.1561 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30664	31879	23.02	49.94
31835	32579	3.07	58.69
31879	31879	0.00	61.19
32579	32579	0.00	58.39
34130	31879	2.85	61.29

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Dextrorotary lactic acid 20%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30740-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Dextrorotary lactic acid 20%; Acidum L (+)-lacticum 20%

Special notes

When selecting the *Dextrorotary lactic acid 20%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dextrorotary lactic acid 20%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Dextrorotary lac...	12105501	30740	40	1405101
Caelo	Dextrorotary lac...	15113303	31720	40	20150622
Caelo	Dextrorotary lac...	162572	33318	40	20170216*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Dextrorotary lactic acid 20%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Dextrorotary lactic acid 20%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dextrorotary lactic acid 20%	15113303	31720 [†]	20
Caelo	Dextrorotary lactic acid 20%	161639	32341	40
Caelo	Dextrorotary lactic acid 20%	172086	33711	40
Caelo	Dextrorotary lactic acid 20%	172086	34083	40

- 93 734 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 4 spectra from 4 *Apo-Ident* customers from 3 batches from the substance/substance group *Dextrorotary lactic acid 20%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Dextrorotary lactic acid 20%	16257201	1
Caelo	Dextrorotary lactic acid 20%	15113303	2
Caelo	Dextrorotary lactic acid 20%	15334302	1

- 91 348 spectra from 501 *Apo-Ident* customers from a total of 27 006 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dextrorotary lactic acid 20%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dextrorotary lactic acid 20%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	140	0	93 734
Type C	0	4	0	91 348

The substance/substance group *Dextrorotary lactic acid 20%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6608 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to

exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30740	30740	0.00	23.34
31720	30740	4.00	21.70
33318	33318	0.00	20.78

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Diethyl ether**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31651-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Diethyl ether

Special notes

When selecting the *Diethyl ether* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Diethyl ether*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Hedinger	Diethyl ether	052159	31854	40	9337
Otto Fischar	Diethyl ether	5975088	31855	40	20150817
Hedinger	Diethyl ether	052162	34122	40	20170216

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Diethyl ether*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 200 spectra from 6 reference samples of the substance/substance group *Diethyl ether*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Otto Fischar	Diethyl ether	5970071	30178	40
Otto Fischar	Diethyl ether	zu erfassen	30681	40
Hedinger	Diethyl ether	052159	31854 [†]	20
Otto Fischar	Diethyl ether	5975088	31855 [†]	20
Hedinger	Diethyl ether	052161	33214	40
Otto Fischar	Diethyl ether	5977054	33633	40

- 93 674 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 136 spectra from 76 *Apo-Ident* customers from 77 batches from the substance/substance group *Diethyl ether*.
- These include spectra of independent samples from 77 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Hedinger	Diethyl ether	52152	1
Fischar	Diethyl ether	5971081	1
Fischar/Gehe	Diethyl ether	5973018	1
Otto Fischar/GEHE	Diethyl ether	597301B/311152/2	1
Fischar/Gehe	Diethyl ether	597405A	1
Hedinger/ Noweda	Diethyl ether	52156	4
Hedinger/Gehe	Diethyl ether	52156	2
Hedinger/Gehe	Diethyl ether	52157	1
Hedinger	Diethyl ether	52156	1
Hedinger/Noweda	Diethyl ether	52157	3
Hedinger/Kehr	Diethyl ether	52157	1
Hedinger	Diethyl ether	52161	1
VWR	Diethyl ether	14B110030	1
Fischar/Noweda	Diethyl ether	5973019	1
Hedinger	Diethyl ether	14070804	1
Otto Fischar/Sanacorp	Diethyl ether	597301A	1
Fischar	Diethyl ether	5973065	1
Fischar /Zentrum-Apotheke	Diethyl ether	5974013	1
Fischar	Diethyl ether	597401A	1
Alliance Healthcare/Fis...	Diethyl ether	5973065	1
Alliance Healthcare/Fis...	Diethyl ether	597401A	1
Otto FischerGmbH / Alli...	Diethyl ether	597401B	1
FISCHER GMBH / GEHE	Diethyl ether	5974051	1
Fischar	Diethyl ether	5974051	1
otto fischar/Sana	Diethyl ether	597401B	1
Caelo	Diethyl ether	5974058	1
noweda/fischar	Diethyl ether	5974056	1
Stumpf/Fischer	Diethyl ether	5974052	1
Fischar	Diethyl ether	5974053	1
Fischar/Phoenix	Diethyl ether	5974058	1
Fischar	Diethyl ether	5974059	2
Otto Fischar GmbH & Co...	Diethyl ether	5974102	2
E+J, 28.09.2016	Diethyl ether	597601C	1
fishar	Diethyl ether	5975012	1
Otto Fischar / Sanacorp	Diethyl ether	5975015	1
Fischer/AEP	Diethyl ether	5976061	1
Otto Fischer/ Kehr	Diethyl ether	5976062	1
Fischar / Gehe	Diethyl ether	5976063	1
Fischar/Sanacorp	Diethyl ether	597606	1
noweda/ Otto Fischar	Diethyl ether	5976066	1
Otto Fischar GmbH & Co...	Diethyl ether	5976065	1
Caelo	Diethyl ether	5976064	1
Allianc. Health. / Otto...	Diethyl ether	597606B	1
Otto Fischar	Diethyl ether	597606B	1
Caelo	Diethyl ether	597606B	1
fishar/ Sanacorp	Diethyl ether	5977054	2
Fischar/Kehr	Diethyl ether	597606B	2
Fischar	Diethyl ether	597606B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischer/Phönix	Diethyl ether	597606b	1
Phönix/Fischer	Diethyl ether	5977018	1
Fischar/Kehr	Diethyl ether	5977055	1
fischar	Diethyl ether	5977056	1
Otto Fischer/Kehr	Diethyl ether	5977056	1
hedinger/anzag	Diethyl ether	52153	2
Otto Fischar / Noweda	Diethyl ether	77016	1
Sanacorp	Diethyl ether	52152	1
Phönix	Diethyl ether	52153	1
Hedinger	Diethyl ether	52153	4
Hedinger/Gehe	Diethyl ether	52153	1
Audor Pharma	Diethyl ether	52153	2
Hedinger/V.d.L.	Diethyl ether	52153	1
Hedinger	Diethyl ether	52155	1
gepepharm GmbH	Diethyl ether	52155	1
Hedinger/Noweda	Diethyl ether	52153	1
Caelo	Diethyl ether	52156	1
Hedinger/Kehr	Diethyl ether	52156	1
Fischar	Diethyl ether	52157	1
Hedinger/Sanacorp	Diethyl ether	12114902	1
Hedinger	Diethyl ether	30011406	1
Caelo	Diethyl ether	28051304	1
Noweda	Diethyl ether	5971381	1
Fischar/Sanacorp	Diethyl ether	5972042	1
Otto Fischar GmbH	Diethyl ether	5972045	1
Fischer Saarbrücken/GEHE	Diethyl ether	5972047	1
Otto Fischar GmbH& Co.KG	Diethyl ether	5972049	1
Otto Fischar GmbH	Diethyl ether	5972048	1
Fischar 001	Diethyl ether	597204B	1
Otto Fischar GmbH& Co.KG	Diethyl ether	597204C	1
Phönix	Diethyl ether	597204B	1
Otto Fischar GmbH	Diethyl ether	597208	1
Gehe	Diethyl ether	5972082	1
Fischar	Diethyl ether	5972082	1
Otto Fischar GmbH	Diethyl ether	5972082	1
fischar	Diethyl ether	5972083	1
Otto Fischar GmbH	Diethyl ether	5972088	1
Otto Fischar GmbH	Diethyl ether	5972089	1
Fischar	Diethyl ether	5972089	1
Fischar	Diethyl ether	597208C	1
Gehe/Fischar	Diethyl ether	597301	2
Otto Fischar GmbH	Diethyl ether	597301B	1
Fischar	Diethyl ether	5973013	1
Caelo	Diethyl ether	5973014	1
Fischar/ Now	Diethyl ether	597301C	1
Phönix/Fischar	Diethyl ether	5973014	1
Fischar/Anzag	Diethyl ether	5973063	1
Fischar/Sanacorp	Diethyl ether	5973066	1
Fischar/Phönix	Diethyl ether	5973064	1
Fischar	Diethyl ether	5973063	1
Phoenix/Fischarl	Diethyl ether	5973066	1
Fischer	Diethyl ether	5973066	1
Fischar	Diethyl ether	5973069	1
Fischar	Diethyl ether	5973067	1
Fischer	Diethyl ether	597306A	1
Fischar	Diethyl ether	5974013	1
Fischar	Diethyl ether	597401C	1
Fischar/Phoenix	Diethyl ether	5974014	1
Otto Fischar	Diethyl ether	5975089	1
fischar	Diethyl ether	597601C	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischar/ Now	Diethyl ether	597306C	1
Hedinger/Anzag	Diethyl ether	52161	1
Hedinger/Noweda	Diethyl ether	52161	2
Fischar/Sanacorp	Diethyl ether	597606B	1
Hedinger	Diethyl ether	52162	1
Fischar/Sanacorp	Diethyl ether	5977018	1
Ichthyol-Gesellschaft	Diethyl ether	5976011	1
Fischar/Sanacorp	Diethyl ether	5977056	3
Fischar/Sanacorp	Diethyl ether	5978014	1

- 91 216 spectra from 501 *Apo-Ident* customers from a total of 26 932 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Diethyl ether* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Diethyl ether* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	190	10	93 674
Type C	0	134	2	91 216

The substance/substance group *Diethyl ether* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	95.0000 % (> 93.5000 %)
Type C	100.0000 % (> 99.6586 %)	98.5294 % (> 96.3235 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31854	31854	0.00	142.36
31855	31854	0.75	142.45
34122	31854	6.05	137.25

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Dimethyl sulfoxide
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31498-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Dimethyl sulfoxide; Dimethylsulfoxidum

Special notes

When selecting the *Dimethyl sulfoxide* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dimethyl sulfoxide*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Euro OTC	Dimethyl sulfoxide	131000802	31498	40	AR-15-FG-002774-01
Fagron	Dimethyl sulfoxide	15A16-C01-16194	31732	40	20150130
Euro OTC	Dimethyl sulfoxide	1508018	31734	40	0000010817
Euro OTC	Dimethyl sulfoxide	1511013	31803	40	20151124*
Euro OTC	Dimethyl sulfoxide	1511028	31813	40	20151202*
Euro OTC	Dimethyl sulfoxide	1601048	31848	40	20160129*
Euro OTC	Dimethyl sulfoxide	1604010	32280	40	20160513*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 280 spectra of 7 reference samples from the substance/substance group *Dimethyl sulfoxide*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 7 different batches.
- 35 264 spectra from a total of 833 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 800 spectra from 23 reference samples of the substance/substance group *Dimethyl sulfoxide*.
- These include spectra of independent samples from 17 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Dimethyl sulfoxide	131000802	31498 [†]	20
Fagron	Dimethyl sulfoxide	15A16-C01-16194	31732 [†]	20
Euro OTC	Dimethyl sulfoxide	1508018	31734 [†]	20
Euro OTC	Dimethyl sulfoxide	1511013	31803 [†]	20
Euro OTC	Dimethyl sulfoxide	1511028	31813 [†]	20
Euro OTC	Dimethyl sulfoxide	1601048	31848 [†]	20
Euro OTC	Dimethyl sulfoxide	1606007	32460	40
Euro OTC	Dimethyl sulfoxide	1608009	32707	40
Euro OTC	Dimethyl sulfoxide	1608029	32708	40
Euro OTC	Dimethyl sulfoxide	1609006	32713	40
Euro OTC	Dimethyl sulfoxide	1610015	32835	40
Euro OTC	Dimethyl sulfoxide	1701006	33233	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Dimethyl sulfoxide	1705025	33451	40
Euro OTC	Dimethyl sulfoxide	1705005	33452	40
Euro OTC	Dimethyl sulfoxide	1709010	33666	40
Euro OTC	Dimethyl sulfoxide	1710001	33667	40
Euro OTC	Dimethyl sulfoxide	1801008	33925	40
Euro OTC	Dimethyl sulfoxide	1711020	33933	40
Euro OTC	Dimethyl sulfoxide	1803015	34097	40
Euro OTC	Dimethyl sulfoxide	1805002	34181	40
Euro OTC	Dimethyl sulfoxide	1806021	34193	40
Euro OTC	Dimethyl sulfoxide	1808012	34216	40
Caelo	Dimethyl sulfoxide	181976	34221	40

- 93 074 spectra from a total of 2144 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 530 spectra from 217 *Apo-Ident* customers from 105 batches from the substance/substance group *Dimethyl sulfoxide*.
- These include spectra of independent samples from 100 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Dimethyl sulfoxide	11H30-N06	3
Fagron	Dimethyl sulfoxide	11K29-N03	3
Hedinger/Noweda	Dimethyl sulfoxide	1006007-01	1
Caelo	Dimethyl sulfoxide	1206010-01	1
Euro OTC	Dimethyl sulfoxide	130008-01	1
Euro OTC	Dimethyl sulfoxide	1306010-01	35
PHOENIX 06.03.2014	Dimethyl sulfoxide	1306010-01	1
Audor Pharma	Dimethyl sulfoxide	1306010-01D	1
Fagron	Dimethyl sulfoxide	1306010-02	1
Caelo	Dimethyl sulfoxide	1306010-02	1
Euro OTC	Dimethyl sulfoxide	1306010-02	22
	Dimethyl sulfoxide	1306010-02	1
Chemical corp. / Noweda	Dimethyl sulfoxide	1306010-02	1
phoenix 06.03.14	Dimethyl sulfoxide	1306010-02	1
Euro OTC	Dimethyl sulfoxide	13063010-01	1
Euro OTC	Dimethyl sulfoxide	131000801	1
Euro OTC	Dimethyl sulfoxide	1310008-01	20
Caelo	Dimethyl sulfoxide	1310008-01	2
Euro OTC	Dimethyl sulfoxide	1310008-02	17
phoenix 06.05.14	Dimethyl sulfoxide	1310008-01	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dimethyl sulfoxide	1310008-02	1
Fagron	Dimethyl sulfoxide	1310008-02	1
Chemical Corp./Sanacorp	Dimethyl sulfoxide	1310008-02	1
EuRho/Phönix	Dimethyl sulfoxide	1310008-02	1
Euro OTC	Dimethyl sulfoxide		3
Phoenix 10.06.14	Dimethyl sulfoxide	131008-01	1
Euro	Dimethyl sulfoxide	140403-01	1
Euro OTC	Dimethyl sulfoxide	1404031-01	10
Caelo	Dimethyl sulfoxide	1404031-01	1
Gaylord/Sanacorp 28.08.14	Dimethyl sulfoxide	1404031-01	1
Noweda	Dimethyl sulfoxide	140403101	1
Chemical Corp.	Dimethyl sulfoxide	1404031-02	1
Phönix Gotha	Dimethyl sulfoxide	1404031-02	1
Euro OTC	Dimethyl sulfoxide	1404031-02	8
Euro/ Gehe	Dimethyl sulfoxide	1404031-02	1
Noweda	Dimethyl sulfoxide	1404031-02	1
Euro OTC	Dimethyl sulfoxide	14051404	1
Gaylord ChemicalCorp./ ...	Dimethyl sulfoxide	1404031-02	1
Euro OTC	Dimethyl sulfoxide	1411039-01	14
Euro OTC	Dimethyl sulfoxide	141039-02	1
Euro OTC	Dimethyl sulfoxide	141103901	1
Euro OTC	Dimethyl sulfoxide	1411039-02	7
	Dimethyl sulfoxide	1411039-02	1
EuRho	Dimethyl sulfoxide	1411039-02	1
Euro OTG/ AHD	Dimethyl sulfoxide	1411039-02	1
Euro OTC	Dimethyl sulfoxide	141103902	1
Caelo	Dimethyl sulfoxide	1411039-02	3
Euro OTC	Dimethyl sulfoxide	14H28-M01	2
Fagron	Dimethyl sulfoxide	14H28-M01	2
Fagron	Dimethyl sulfoxide	14H28-M02	1
Fagron	Dimethyl sulfoxide	14H28-M03	1
Euro OTC	Dimethyl sulfoxide	1505021	1
Euro OTC	Dimethyl sulfoxide	15011503	1
Phönix Gotha	Dimethyl sulfoxide	1505021-01	1
Euro OTC	Dimethyl sulfoxide	1505021-01	17
Gaylord / Fiebig	Dimethyl sulfoxide	1505021-01	1
Euro OTC	Dimethyl sulfoxide	1508013-01	1
Euro OTC	Dimethyl sulfoxide	1508018-01	22
Gaylord Chemical Corp./...	Dimethyl sulfoxide	1508018-01	1
Caelo	Dimethyl sulfoxide	1508018-01	2
Caerol / Gehe	Dimethyl sulfoxide	1508018-01	1
Euro OTC	Dimethyl sulfoxide	150801801	1
Sanacorp	Dimethyl sulfoxide	1508018-01	1
Euro OTC	Dimethyl sulfoxide	1508018-02	20
Caelo	Dimethyl sulfoxide	1508018-02	1
Gaylord Chemical Corp./...	Dimethyl sulfoxide	1508018-02	2
Euro OTC	Dimethyl sulfoxide	1511028-01	3
AHD	Dimethyl sulfoxide	151102801	1
Euro OTC	Dimethyl sulfoxide	151102801	2
Fagron	Dimethyl sulfoxide	15A16-C01-16194	2
Phönix	Dimethyl sulfoxide	1511028-01	1
Euro OTC	Dimethyl sulfoxide	15111028-01	1
Euro OTC	Dimethyl sulfoxide	160048-02	1
Fagron	Dimethyl sulfoxide	15A16-C01-16195	1
vitalundfitmit100/ Sana...	Dimethyl sulfoxide		1
Euro OTC	Dimethyl sulfoxide	1601048-01	8
Phoenix	Dimethyl sulfoxide	1601048-01	1
Caelo	Dimethyl sulfoxide	1601048-02	1
Bombastus	Dimethyl sulfoxide	1601048-02	1
Euro OTC	Dimethyl sulfoxide	1601048-02	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alliance	Dimethyl sulfoxide	1601048-02	1
Euro OTC	Dimethyl sulfoxide	1604010-01	11
Ebert 09.08.16 KKH	Dimethyl sulfoxide	1604010-01	1
Ebert 12.08.16	Dimethyl sulfoxide	1604040-01	1
Euro OTC	Dimethyl sulfoxide	1606007-01	13
phönix	Dimethyl sulfoxide	1606007-01	1
Alliance	Dimethyl sulfoxide	1606007-01	1
Alte Apotheke	Dimethyl sulfoxide	1606007-01	1
Sanacorp	Dimethyl sulfoxide	160600701	1
Ebert KKH29.09.16	Dimethyl sulfoxide	16060074-01	1
Euro OTC	Dimethyl sulfoxide	1609006-01	11
Caelo	Dimethyl sulfoxide	1609006-01	2
Euro OTC	Dimethyl sulfoxide	1610015-01	14
Noweda	Dimethyl sulfoxide	1610015-01	2
Euro OTC	Dimethyl sulfoxide	1610015-02	25
GEHE	Dimethyl sulfoxide	1610015-01	1
EuRho/Noweda	Dimethyl sulfoxide	1610015-02	1
Caelo	Dimethyl sulfoxide	1610015-02	1
Gehe	Dimethyl sulfoxide	1610015-1	1
Euro OTC	Dimethyl sulfoxide	161005-01	1
Euro OTC	Dimethyl sulfoxide	1705025-01	8
Euro OTG/ AHD	Dimethyl sulfoxide	1705025-01	3
Noweda	Dimethyl sulfoxide	170502501	1
Fagron	Dimethyl sulfoxide	1705025-01	1
Euro OTC	Dimethyl sulfoxide	1709010-01	2
Caelo	Dimethyl sulfoxide	170901001	1
Euro OTC	Dimethyl sulfoxide	170522DS/1610015-01	1
Euro OTC	Dimethyl sulfoxide	171000-01	1
Sanacorp W.E: 23.01.18 ...	Dimethyl sulfoxide	1710001-01	1
Euro OTC	Dimethyl sulfoxide	1711020-01	3
Euro OTC	Dimethyl sulfoxide	1805002-01	1
Euro OTG/ AHD	Dimethyl sulfoxide	1711020-01	1
Audor Pharma	Dimethyl sulfoxide	AP1306010-01D	3
Euro OTC	Dimethyl sulfoxide	1813A-01336	1
Euro OTC	Dimethyl sulfoxide	SYN140623	1
Euro OTC	Dimethyl sulfoxide	SYN160215	1
EuRho	Dimethyl sulfoxide	SYN160216	1
Euro OTC	Dimethyl sulfoxide	2071301	1
Euro OTC	Dimethyl sulfoxide	2071302	1
Euro OTC	Dimethyl sulfoxide	5111303	1
Euro OTC	Dimethyl sulfoxide	613A-01336	1
Euro OTC	Dimethyl sulfoxide	1004041-01	1
Alliance Health 06.09.2...	Dimethyl sulfoxide	10G26-N12	1
Euro OTC	Dimethyl sulfoxide	1106006-01	4
Euro OTC	Dimethyl sulfoxide	110600601	1
Euro OTC	Dimethyl sulfoxide	1111059-01	1
ANZAG am 21.05.2012	Dimethyl sulfoxide	1111059-01	2
Noweda	Dimethyl sulfoxide	1111059-01	1
Gehe	Dimethyl sulfoxide	1111059-01	1
Audor Pharma	Dimethyl sulfoxide	1111059-01D	1
Gehe	Dimethyl sulfoxide	1111059-02	1
Noweda	Dimethyl sulfoxide	1111059-02	1
Phönix Gotha	Dimethyl sulfoxide	1111059-02	1
EuRho	Dimethyl sulfoxide	1111059-02	1
Fagron	Dimethyl sulfoxide	11E03-M04	2
Fagron	Dimethyl sulfoxide	11A21-N03	1
Euro OTC	Dimethyl sulfoxide	1111059-02	1
Fagron	Dimethyl sulfoxide	11G18-N15	1
Euro OTC	Dimethyl sulfoxide	1207008-01	2
Sanacorp	Dimethyl sulfoxide	1207008-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Dimethyl sulfoxide	1207008-02	7
Euro OTC	Dimethyl sulfoxide	1210059-01	21
Großhandel	Dimethyl sulfoxide	1210059-01	1
EuRho	Dimethyl sulfoxide	1210059-01	1
Alliance Healthcare	Dimethyl sulfoxide	1210059-01	1
gaylord	Dimethyl sulfoxide	1210059-01	1
Euro OTC	Dimethyl sulfoxide	121005901	1
Gehe 18.06.2013	Dimethyl sulfoxide	1210059-01	1
Phönix	Dimethyl sulfoxide	1210059-01	1
Euro OTC	Dimethyl sulfoxide	12101202	1
Phoenix 23.08.2012	Dimethyl sulfoxide	12a09-N03	1
Euro OTC	Dimethyl sulfoxide	13003010-01	1
Euro OTC	Dimethyl sulfoxide	13011403	1
Euro OTC	Dimethyl sulfoxide	1303010-01	17
Caelo	Dimethyl sulfoxide	1306010-01	2
Gaylord/Sanacorp	Dimethyl sulfoxide	1306010-01	1
Alliance Health 04.11.2...	Dimethyl sulfoxide	1306010-01	1
AHD,17.10.13	Dimethyl sulfoxide	1306010-02	1
Euro/ Fiebig	Dimethyl sulfoxide	1306010-02	1
gehe 07.01.2014	Dimethyl sulfoxide	1306010-02	2
Fagron	Dimethyl sulfoxide	13071201	1
Caeöo	Dimethyl sulfoxide	1310008-01	1
Phönix	Dimethyl sulfoxide	1310008-01	1
Caelo	Dimethyl sulfoxide	131000801	1
Euro OTC	Dimethyl sulfoxide	1404031	1
Fagron	Dimethyl sulfoxide	1404031-02	1
Euro OTC	Dimethyl sulfoxide	26011301	1
EuRho/Noweda	Dimethyl sulfoxide	1601048-01	2
Fagron	Dimethyl sulfoxide	27061203	1
Caelo	Dimethyl sulfoxide	29081202	1
Euro OTC	Dimethyl sulfoxide	PHS111205	1
Audor Pharma	Dimethyl sulfoxide	AP1303010-01D	1
Audor Pharma	Dimethyl sulfoxide	ap130301001d	1
Phönix	Dimethyl sulfoxide	AP1303010-01D	1
EuRho/Fiebig	Dimethyl sulfoxide	SYN150525	1
Fagron	Dimethyl sulfoxide	PHS1000510	1
Sanacorp/ PZN08597061	Dimethyl sulfoxide	1306010-01-339413100	1
Sanacorp	Dimethyl sulfoxide	1306010-02	1
Audor Pharma	Dimethyl sulfoxide	AP1210059-01D	1
Phoenix WE 05.04.17 -7,...	Dimethyl sulfoxide	1610015-01	2
Euro OTC	Dimethyl sulfoxide	1606007	1
Phönix	Dimethyl sulfoxide	1705025-01	1
Phönix/EuRho	Dimethyl sulfoxide	1207008-01	1

- 90 822 spectra from 501 *Apo-Ident* customers from a total of 26 905 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dimethyl sulfoxide* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dimethyl sulfoxide* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The

following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	280	0	35 264
Type B	0	780	20	93 074
Type C	0	508	22	90 822

The substance/substance group *Dimethyl sulfoxide* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.8571 %)
Type B	100.0000 % (> 99.9761 %)	97.5000 % (> 97.1250 %)
Type C	100.0000 % (> 99.6586 %)	95.8491 % (> 95.2830 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31498	31498	0.00	146.91
31732	31498	1.61	150.80
31734	31734	0.00	150.83
31803	31803	0.00	151.94
31813	31813	0.00	151.87
31848	31848	0.00	150.85
32280	32280	0.00	155.96

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Dimeticon 100 / 350 / 1000**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30660-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Dimeticon 100 / 350 / 1000; Dimeticone 100; Dimeticone 1000; Dimeticone 350; Dimeticonum 100; Dimeticonum 1000; Dimeticonum 350

Special notes

When selecting the *Dimeticon 100 / 350 / 1000* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dimeticon 100 / 350 / 1000*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Dimeticone 100	12140101	30660	40	1403158
Caelo	Dimeticone 100	13049608	31333	40	20130301
Caelo	Dimeticone 350	13236306	31334	40	AR-15-FG-009277-01
Fagron	Dimeticone 350	14C20-B05-292285	31659	40	20140331

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Dimeticone 1000	161307	32169	40	20160427*
Caelo	Dimeticone 1000	15349907	32587	40	20151106
Euro OTC	Dimeticone 350	1702002	33234	40	20170302*
Fagron	Dimeticone 350	14D30-B07-330851	33278	40	20140513
Caelo	Dimeticone 1000	16348504	33654	40	20170111
Caelo	Dimeticone 100	172164	33918	40	20170809*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 400 spectra of 10 reference samples from the substance/substance group *Dimeticone 100 / 350 / 1000*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 10 different batches.
- 35 144 spectra from a total of 830 batches from further 209 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 280 spectra from 8 reference samples of the substance/substance group *Dimeticone 100 / 350 / 1000*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dimeticone 350	10111314	30245	40
Fagron	Dimeticone 350	14C20-B05-292285	31659 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dimeticone 1000	161307	32169 [†]	20
Caelo	Dimeticone 350	160754	32208	40
Caelo	Dimeticone 100	15199409	32304	40
Fagron	Dimeticone 1000	14C20-B04-306781	32485	40
Caelo	Dimeticone 100	15404609	33628	40
Caelo	Dimeticone 1000	171981	33919	40

- 93 594 spectra from a total of 2159 batches from further 434 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 113 spectra from 53 *Apo-Ident* customers from 53 batches from the substance/substance group *Dimeticone 100 / 350 / 1000*.
- These include spectra of independent samples from 48 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Dimeticone 350	13236306	13
Caelo	Dimeticone 350	13236302	2
Caelo	Dimeticone 350	13049802	2
Caelo	Dimeticone 350	14074509	15
Euro OTC	Dimeticone 350	1407004-01	1
Jenne	Dimeticone 350	14074509	1
Caelo	Dimeticone 350	14074512	9
Caelo	Dimeticone 350	14704509	1
Caelo	Dimeticone 350	14074518	1
Caelo	Dimeticone 350	14336402	2
Ph.Eur.8.0	Dimeticone 350	14C20B05292285	1
Fagron	Dimeticone 350	14C20-B05-292285	5
Euro OTC	Dimeticone 350	14C20-B05-292285	1
Fagron	Dimeticone 350	14D30-B07-300432	1
Caelo	Dimeticone 350	14D30-B07-300432	1
Caelo	Dimeticone 350	15349907	1
Fagron	Dimeticone 350	16075401	1
Caelo	Dimeticone 350	16075407	2
Caelo	Dimeticone 1000	16348503	1
AHD	Dimeticone 350	17090204	1
EuRo/Phönix	Dimeticone 350	1207015-01	1
Beiersdorf Ag	Dimeticone 350	1207015-02	1
Caelo	Dimeticone 350	13236308	2
Caelo	Dimeticone 350	14336403	1
Caelo	Dimeticone 350	13049803	2

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dimeticone 350	13236307	1
Fargon	Dimeticone 350	14C20-B05292282	1
Caelo	Dimeticone 100	15199401	1
Caelo	Dimeticone 350	16075401	1
Caelo	Dimeticone 350	16075405	1
Caelo	Dimeticone 350	17090201	1
Caelo	Dimeticone 350	13049813	1
Caelo	Dimeticone 350	13236312	2
Gehe	Dimeticone 350	1407004-01	1
Caelo	Dimeticone 100	13362904	1
Caelo	Dimeticone 350	14074511	1
Caelo	Dimeticone 100	11215905	1
Caelo	Dimeticone 100	12140101	1
Fagron	Dimeticone 350	14C20-B05	1
Caelo	Dimeticone 100	12310801	1
Caelo	Dimeticone 100	13049603	1
Noweda	Dimeticone 100	12316701	1
Caelo	Dimeticone 350	9371810	1
Fagron	Dimeticone 350	9B10	1
Caelo	Dimeticone 350	10111314	1
Caelo	Dimeticone 350	11042104	2
Großhandel	Dimeticone 350	11317804	1
Kehr	Dimeticone 350	1207015-01	1
Euro OTC	Dimeticone 350	1207015-02	1
Euro OTC	Dimeticone 350	1207015-01	1
Caelo	Dimeticone 350	12186301	1
Anzag 24.01.2013	Dimeticone 350	12186301	1
Caelo	Dimeticone 350	12193302	1
Caelo	Dimeticone 350	12193303	2
Caelo	Dimeticone 350	12193307	2
Caelo	Dimeticone 350	12193313	1
Caelo	Dimeticone 350	12193308	1
Caelo	Dimeticone 350	12193315	1
Caelo	Dimeticone 350	13049804	3
AHD, 07.11.13	Dimeticone 350	13049804	1
Caelo	Dimeticone 350	13049809	1
Fagron	Dimeticone 350	3111Q-02251	1
Caelo	Dimeticone 350		1

- 91 239 spectra from 501 *Apo-Ident* customers from a total of 26 957 batches from a further 331 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dimeticon 100 / 350 / 1000* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dimeticon 100 / 350 / 1000* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	400	0	35 144
Type B	0	280	0	93 594
Type C	0	109	4	91 239

The substance/substance group *Dimeticon 100 / 350 / 1000* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.8571 %)
Type C	100.0000 % (> 99.6586 %)	96.4602 % (> 93.8053 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30660	30660	0.00	431.85
31333	32169	1.58	439.64
31334	31334	0.00	441.82
31659	33918	2.32	438.22
32169	32169	0.00	449.44
32587	33234	2.72	428.94
33234	33234	0.00	428.04
33278	33234	1.94	424.79
33654	33918	2.24	436.79
33918	33918	0.00	424.84

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Dwarf pine oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30204-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Dwarf pine oil; Oleum pini pumilionis; Pini pumilionis aetherolum

Special notes

When selecting the *Dwarf pine oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dwarf pine oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Dwarf pine oil	12C29-N17	31184	40	1403167
Caelo	Dwarf pine oil	11375614	31189	40	20120102
Euro OTC	Dwarf pine oil	1210024-02	31529	40	20140806

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Euro OTC	Dwarf pine oil	1511014	31818	40	20151214*
Caelo	Dwarf pine oil	16114705	33266	40	20160420
Caelo	Dwarf pine oil	16114706	33282	40	20160420

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 240 spectra of 6 reference samples from the substance/substance group *Dwarf pine oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 6 different batches.
- 35 304 spectra from a total of 834 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 320 spectra from 9 reference samples of the substance/substance group *Dwarf pine oil*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dwarf pine oil	10060307	30204	40
Euro OTC	Dwarf pine oil	1210024-02	31529 [†]	20
Euro OTC	Dwarf pine oil	1511014	31818 [†]	20
Euro OTC	Dwarf pine oil	1606003	32470	40
Euro OTC	Dwarf pine oil	1703031	33235	40
Fagron	Dwarf pine oil	15A05-B06-322743	33276	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dwarf pine oil	172158	33679	40
Euro OTC	Dwarf pine oil	1801001	33931	40
Euro OTC	Dwarf pine oil	1807002	34194	40

- 93 554 spectra from a total of 2158 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 160 spectra from 80 *Apo-Ident* customers from 71 batches from the substance/substance group *Dwarf pine oil*.
- These include spectra of independent samples from 66 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Dwarf pine oil	621H23	1
Phönix	Dwarf pine oil	11/19-N12	1
Caelo	Dwarf pine oil	11375614	3
Taoasis	Dwarf pine oil	1148-148	1
Euro OTC	Dwarf pine oil	1210024-02	1
Frey+Lau GmbH/Phönix	Dwarf pine oil	12C29-N17	2
Fagron	Dwarf pine oil	12-K06-N06	1
Fagron	Dwarf pine oil	12K06-N06	1
Caelo	Dwarf pine oil	13084309	1
Fagron	Dwarf pine oil	13C04-N02	7
Frey & Lau GmbH/Phoenix	Dwarf pine oil	13C04-N020133/265	1
Fagron	Dwarf pine oil	13E06-N04	5
Fagron	Dwarf pine oil	13H28-N04	3
Caelo	Dwarf pine oil	14083301	7
Caelo	Dwarf pine oil	14083302	3
Caelo	Dwarf pine oil	14083303	1
Caelo	Dwarf pine oil	14083311	1
Caelo	Dwarf pine oil	14083315	1
Caelo	Dwarf pine oil	14083310	1
Euro OTC	Dwarf pine oil	1411021-01	5
Fagron	Dwarf pine oil	15051509	1
Euro OTC	Dwarf pine oil	1511014-01	2
Caelo	Dwarf pine oil	15239201	3
Caelo	Dwarf pine oil	15239202	2
Caelo	Dwarf pine oil	153239202	1
Caelo	Dwarf pine oil	15239203	1
Caelo	Dwarf pine oil	15375501	6
Caelo	Dwarf pine oil	15375502	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dwarf pine oil	15375505	4
Fagron	Dwarf pine oil	15A05B06	1
Fagron	Dwarf pine oil	15A05-B06	5
Fagron	Dwarf pine oil	15A05-B06-305064	1
Euro OTC	Dwarf pine oil	1606003-01	1
Caelo	Dwarf pine oil	16114701	5
Caelo	Dwarf pine oil	1606003-01	1
Ichthyol-Gesellschaft	Dwarf pine oil	1610280P/15375501	1
Fagron	Dwarf pine oil	15A05-B06-305065	1
Caelo	Dwarf pine oil	16114705	6
Caelo	Dwarf pine oil	17215801	2
Caelo	Dwarf pine oil	17215802	1
Taoasis	Dwarf pine oil	25599-117839	1
Caelo	Dwarf pine oil	21051403	1
Taoasis	Dwarf pine oil	25599-118829	1
Fritz Schiele	Dwarf pine oil		1
Taoasis	Dwarf pine oil	586/116867	1
Fagron	Dwarf pine oil	736339	3
Taoasis	Dwarf pine oil	586-116867	1
Taoasis	Dwarf pine oil	8783-125163	1
Caelo	Dwarf pine oil	11069807	1
Caelo	Dwarf pine oil	10020505	1
Caelo	Dwarf pine oil	1792977	2
Caelo	Dwarf pine oil	11357601	1
Caelo	Dwarf pine oil	11375601	11
Klenk	Dwarf pine oil	11375601	1
Alliance	Dwarf pine oil	11375601	1
Caelo	Dwarf pine oil	11375608	1
Caelo	Dwarf pine oil	11375613	3
Fagron	Dwarf pine oil	11G18-N10	2
Fagron	Dwarf pine oil	11L19-N12	1
Caelo	Dwarf pine oil	12239902	2
Caelo	Dwarf pine oil	12239904	3
Caelo	Dwarf pine oil	12239903	2
Fagron	Dwarf pine oil	12C29-N17	2
Fagron	Dwarf pine oil	12G12-N11	2
Fargo	Dwarf pine oil	13C04-N02	1
PVL	Dwarf pine oil	1360313	1
Fagron	Dwarf pine oil	13C04-N02	1
Caelo	Dwarf pine oil	17012014	2
Caelo	Dwarf pine oil	23021301	1
Fagron	Dwarf pine oil	27111207	1
Taoasis GmbH	Dwarf pine oil	30985-110461	1
Taoasis	Dwarf pine oil	30985-110461	1
Taoasis	Dwarf pine oil	30985-111105	2
Taoasis	Dwarf pine oil	33347-112503	1
Fagron	Dwarf pine oil	14D14-B01	6
VDL; 7.50EUR; 28.9.2015	Dwarf pine oil	15239201	1
Caelo	Dwarf pine oil	1137575601	1
Taoasis	Dwarf pine oil	30985-106654	1

- 91 192 spectra from 501 *Apo-Ident* customers from a total of 26 940 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dwarf pine oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dwarf pine oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	240	0	35 304
Type B	0	320	0	93 554
Type C	0	150	10	91 192

The substance/substance group *Dwarf pine oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.1250 %)
Type C	100.0000 % (> 99.6586 %)	93.7500 % (> 91.8750 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31184	31184	0.00	30.91
31189	31184	4.44	31.99
31529	31818	2.12	25.98
31818	31818	0.00	26.47
33266	31818	6.18	28.79
33282	31818	6.93	30.36

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Ethanol 90% pure / denatured**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30379-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Ethanol 90% pure / denatured; Ethanolum 90% rein / vergällt

Special notes

When selecting the *Ethanol 90% pure / denatured* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Ethanol 90% pure / denatured*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
apomix	Ethanol 90% pure...	25A040913	31056	40	1404654
Caelo	Ethanol 90% pure...	14052802	31392	40	20140225
Caelo	Ethanol 90% pure...	163144	33102	40	20161209*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Ethanol 90% pure / denatured*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 280 spectra from 7 reference samples of the substance/substance group *Ethanol 90% pure / denatured*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Ethanol 90% pure / denatured	11069005	30319	40
Caelo	Ethanol 90% pure / denatured	11329301	30379	40
Caelo	Ethanol 90% pure / denatured	161944	32324	40
Caelo	Ethanol 90% pure / denatured	163142	33103	40
Caelo	Ethanol 90% pure / denatured	163143	33104	40
Caelo	Ethanol 90% pure / denatured	171741	33714	40
Caelo	Ethanol 90% pure / denatured	171739	33715	40

- 93 594 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 461 spectra from 165 *Apo-Ident* customers from 264 batches from the substance/substance group *Ethanol 90% pure / denatured*.
- These include spectra of independent samples from 259 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Zentrum-Apotheke	Ethanol 90% pure / denatured	514I-03042	1
Caelo	Ethanol 90% pure / denatured	5A130214	1
Apomix	Ethanol 90% pure / denatured	2A130114	1
Apotärze/Zentrum	Ethanol 90% pure / denatured	10092013-D	1
Caelo	Ethanol 90% pure / denatured	814E-03042	1
Anzag am 29.01.2013	Ethanol 90% pure / denatured	12068403	2
Caelo	Ethanol 90% pure / denatured	12275402	2
Caelo	Ethanol 90% pure / denatured	12355401	4
Anzag 28.03.2013	Ethanol 90% pure / denatured	12355601	1
Caelo	Ethanol 90% pure / denatured	12355501	2
Caelo	Ethanol 90% pure / denatured	13019602	2
Caelo	Ethanol 90% pure / denatured	13019601	1
Caelo	Ethanol 90% pure / denatured	12355601	4
Caelo	Ethanol 90% pure / denatured	13019701	2
Phoenix	Ethanol 90% pure / denatured	13019703	1
Caelo	Ethanol 90% pure / denatured	13019703	2
Caelo	Ethanol 90% pure / denatured	13085901	4
Caelo	Ethanol 90% pure / denatured	13086002	1
Caelo	Ethanol 90% pure / denatured	13086001	4
Caelo	Ethanol 90% pure / denatured	13160501	2
Caelo	Ethanol 90% pure / denatured	1305801	1
Caelo	Ethanol 90% pure / denatured	13160601	4
Caelo	Ethanol 90% pure / denatured	13160701	2
Caelo	Ethanol 90% pure / denatured	13290201	4
Caelo	Ethanol 90% pure / denatured	13160703	2
Caelo	Ethanol 90% pure / denatured	13290202	2
Caesar & Loretz GmbH, 4...	Ethanol 90% pure / denatured	13290202	1
Noweda	Ethanol 90% pure / denatured	13290203	1
Caelo	Ethanol 90% pure / denatured	13290203	1
Caelo	Ethanol 90% pure / denatured	13290402	3
Phönix	Ethanol 90% pure / denatured	13399604	1
Caelo	Ethanol 90% pure / denatured	13290301	2
Caelo	Ethanol 90% pure / denatured	13399601	1
Caelo	Ethanol 90% pure / denatured	13399606	2
Caelo	Ethanol 90% pure / denatured	13399605	4
Caelo	Ethanol 90% pure / denatured	13399801	1
Caelo	Ethanol 90% pure / denatured	14027502	1
Caelo	Ethanol 90% pure / denatured	14208901	3
Phonix	Ethanol 90% pure / denatured	14052601	1
Post-Apotheke 15.01.2014	Ethanol 90% pure / denatured	15012014	1
Caelo	Ethanol 90% pure / denatured	15010602	1
Phoenix	Ethanol 90% pure / denatured	13399802	1
Caelo	Ethanol 90% pure / denatured	14208702	2
Caelo	Ethanol 90% pure / denatured	15219501	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ethanol 90% pure / denatured	15442801	3
Caelo	Ethanol 90% pure / denatured	15442502	2
Caelo	Ethanol 90% pure / denatured	1616902	1
Post-Apotheke	Ethanol 90% pure / denatured	1713Q-03042	1
Gütte/Phönix	Ethanol 90% pure / denatured	20112013-M	1
Caelo	Ethanol 90% pure / denatured	21105102	1
Apomix	Ethanol 90% pure / denatured	37A101212	1
Apotärze	Ethanol 90% pure / denatured	4513M-02303	1
Hirsch	Ethanol 90% pure / denatured	5013Q-02303	3
Martinus-Apotheke	Ethanol 90% pure / denatured	A-160913-2D	2
phönix	Ethanol 90% pure / denatured		2
Caelo	Ethanol 90% pure / denatured	AKD6	2
Martinus-Apotheke	Ethanol 90% pure / denatured	E-080114-1D	2
Hustadt Apotheke	Ethanol 90% pure / denatured	Defektur	1
Eigenherstellung	Ethanol 90% pure / denatured	4161112	1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-040414-1D	2
Mühlen-Apotheke	Ethanol 90% pure / denatured	213Q-02303	1
Herstellung	Ethanol 90% pure / denatured	110416P1	1
Caelo	Ethanol 90% pure / denatured	11329304	1
Caelo	Ethanol 90% pure / denatured	11329402	1
Gehe	Ethanol 90% pure / denatured	12068201	1
Caesar & Loretz GmbH, 4...	Ethanol 90% pure / denatured	12068301	1
Caelo	Ethanol 90% pure / denatured	12144401	1
Caelo	Ethanol 90% pure / denatured	12068301	1
Caelo	Ethanol 90% pure / denatured	12275301	1
Caesar & Loretz GmbH, 4...	Ethanol 90% pure / denatured	12275402	1
Phönix 27.08.2012	Ethanol 90% pure / denatured	12144402	1
Caelo	Ethanol 90% pure / denatured	12355504	1
Phönix Gotha	Ethanol 90% pure / denatured	20A030712	1
Anzag 06.10.12	Ethanol 90% pure / denatured	25A290812	1
Fagron	Ethanol 90% pure / denatured	270614P1	1
Mühlen-Apotheke	Ethanol 90% pure / denatured	4012M-02303	1
apomix	Ethanol 90% pure / denatured	22a080911	1
Apomix	Ethanol 90% pure / denatured	32A141211	1
Caelo	Ethanol 90% pure / denatured	12068202	1
Gutenberg-Apotheke	Ethanol 90% pure / denatured	Defektur	1
Caelo	Ethanol 90% pure / denatured	12275403	1
Apomix	Ethanol 90% pure / denatured	35A211112	1
Caelo	Ethanol 90% pure / denatured	10250602	1
Caelo	Ethanol 90% pure / denatured	13085801	1
Caelo	Ethanol 90% pure / denatured	13290401	2
Caelo	Ethanol 90% pure / denatured	14208801	5
Apomix	Ethanol 90% pure / denatured	21A310713	1
Eigenherstellung	Ethanol 90% pure / denatured	4733	1
Caelo	Ethanol 90% pure / denatured	15150701	6
Caelo	Ethanol 90% pure / denatured	14027401	4
Eigenherstellung	Ethanol 90% pure / denatured	217I-02303	1
Eigenherstellung	Ethanol 90% pure / denatured	115Q-07057	1
Apomix / Phönix	Ethanol 90% pure / denatured	5A130214	1
Apomix / Phoenix	Ethanol 90% pure / denatured	5A150216	1
Caelo	Ethanol 90% pure / denatured	713E-02303	1
apomix	Ethanol 90% pure / denatured	7A170215	1
apomix	Ethanol 90% pure / denatured	8A170215	1
Eigen	Ethanol 90% pure / denatured	817Q-03042	1
Apomix/Kehr	Ethanol 90% pure / denatured	8A200317	1
Hirsch	Ethanol 90% pure / denatured	11233	1
PKH Halle	Ethanol 90% pure / denatured	9A270416	1
HirschApotheke	Ethanol 90% pure / denatured	10151-02303	1
Noweda	Ethanol 90% pure / denatured	10A170314	1
Brinn- Holtz	Ethanol 90% pure / denatured	1258-06/17	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Kastner (eigen)	Ethanol 90% pure / denatured	1215I-02303	1
Apomix	Ethanol 90% pure / denatured	12A150517	1
Kastner	Ethanol 90% pure / denatured	1314I-03042	1
Caelo	Ethanol 90% pure / denatured	13399701	1
Caelo	Ethanol 90% pure / denatured	13399802	1
Caelo	Ethanol 90% pure / denatured	13399803	2
Caelo	Ethanol 90% pure / denatured	14027301	1
Apomix / Alliance Healt...	Ethanol 90% pure / denatured	13A170516	1
Caelo	Ethanol 90% pure / denatured	14027302	1
Caelo	Ethanol 90% pure / denatured	14027303	3
Caelo	Ethanol 90% pure / denatured	14027402	1
Caelo	Ethanol 90% pure / denatured	14027501	2
GEHE 29.04.2014	Ethanol 90% pure / denatured	14027502	2
Caelo	Ethanol 90% pure / denatured	14052601	3
Caelo	Ethanol 90% pure / denatured	14052701	2
Caelo	Ethanol 90% pure / denatured	14052801	3
Caelo	Ethanol 90% pure / denatured	14052802	2
Caelo	Ethanol 90% pure / denatured	14060201	5
Caelo	Ethanol 90% pure / denatured	14060301	3
Caelo	Ethanol 90% pure / denatured	14060302	1
Caelo	Ethanol 90% pure / denatured	14122101	2
Sanacorp 24.11.14 EK 63...	Ethanol 90% pure / denatured	14122101	1
Caelo	Ethanol 90% pure / denatured	14122201	4
Caelo	Ethanol 90% pure / denatured	14122202	1
Caelo	Ethanol 90% pure / denatured	14122302	4
Caelo	Ethanol 90% pure / denatured	14122303	1
Caelo	Ethanol 90% pure / denatured	14122401	3
Caelo	Ethanol 90% pure / denatured	14208802	3
Caelo	Ethanol 90% pure / denatured	15010501	3
Apomix	Ethanol 90% pure / denatured	14A120514	1
Caelo	Ethanol 90% pure / denatured	15010504	1
Caelo	Ethanol 90% pure / denatured	15010502	1
Caelo	Ethanol 90% pure / denatured	15010503	1
Caelo	Ethanol 90% pure / denatured	15010601	1
Caelo	Ethanol 90% pure / denatured	15014001	1
Caelo	Ethanol 90% pure / denatured	15014002	7
Stern	Ethanol 90% pure / denatured	150622004	1
Römer-Apotheke	Ethanol 90% pure / denatured	150428003	1
Caelo	Ethanol 90% pure / denatured	15073202	1
Caelo	Ethanol 90% pure / denatured	15073201	4
Caelo	Ethanol 90% pure / denatured	15073203	2
Caelo	Ethanol 90% pure / denatured	15080301	3
Caelo	Ethanol 90% pure / denatured	15150704	3
Caelo	Ethanol 90% pure / denatured	15150705	1
Libellen Apotheke	Ethanol 90% pure / denatured	1515Q-02303	1
Caelo	Ethanol 90% pure / denatured	15219302	2
Eigenherstellung	Ethanol 90% pure / denatured	1515I-02303	1
Noweda	Ethanol 90% pure / denatured	15219302	1
Caelo	Ethanol 90% pure / denatured	15328201	6
Caelo	Ethanol 90% pure / denatured	15328301	12
Caelo	Ethanol 90% pure / denatured	15328202	2
Noweda	Ethanol 90% pure / denatured	15328301	1
Caelo	Ethanol 90% pure / denatured	15328401	2
Caelo	Ethanol 90% pure / denatured	15442601	4
Caelo	Ethanol 90% pure / denatured	15442602	3
Noweda	Ethanol 90% pure / denatured	15442602	1
Fagron	Ethanol 90% pure / denatured	15442602	1
Caelo	Ethanol 90% pure / denatured	15442702	5
Caelo	Ethanol 90% pure / denatured	15442603	1
Noweda	Ethanol 90% pure / denatured	15442603	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ethanol 90% pure / denatured	15442701	1
Caelo	Ethanol 90% pure / denatured	15442703	2
Caelo	Ethanol 90% pure / denatured	16106802	3
Caelo	Ethanol 90% pure / denatured	16106803	3
Caelo	Ethanol 90% pure / denatured	16107003	3
Audor Pharma	Ethanol 90% pure / denatured	1614E-03042	1
Gehe	Ethanol 90% pure / denatured	16106901	1
Caelo	Ethanol 90% pure / denatured	16194402	1
Caelo	Ethanol 90% pure / denatured	16194401	1
Caelo	Ethanol 90% pure / denatured	16194403	3
Caelo	Ethanol 90% pure / denatured	16194501	2
Sanacorp	Ethanol 90% pure / denatured	16194601	1
Caelo	Ethanol 90% pure / denatured	16194601	3
Caelo	Ethanol 90% pure / denatured	16194701	4
Sanacorp	Ethanol 90% pure / denatured	16194501	1
Caelo	Ethanol 90% pure / denatured	16194702	3
Caelo	Ethanol 90% pure / denatured	16314201	4
Caelo	Ethanol 90% pure / denatured	16314302	3
Caelo	Ethanol 90% pure / denatured	16314202	1
AHD	Ethanol 90% pure / denatured	16314401	1
Caelo	Ethanol 90% pure / denatured	16314401	1
Phönix	Ethanol 90% pure / denatured	1694702	1
Apomix/Phönix	Ethanol 90% pure / denatured	16A210610	1
Caelo	Ethanol 90% pure / denatured	17047802	3
Caelo	Ethanol 90% pure / denatured	17047804	2
Caelo	Ethanol 90% pure / denatured	17047901	3
Caelo	Ethanol 90% pure / denatured	17048001	2
Caelo	Ethanol 90% pure / denatured	17048004	2
Eigenherstellung	Ethanol 90% pure / denatured	1716Q-02303	1
DEFEKTUR. Schwanen-Apoth...	Ethanol 90% pure / denatured	170620-1	1
Caelo	Ethanol 90% pure / denatured	17174002	1
Caelo	Ethanol 90% pure / denatured	17173902	1
Caelo	Ethanol 90% pure / denatured	17174102	1
Caelo	Ethanol 90% pure / denatured	17206401	3
Caelo	Ethanol 90% pure / denatured	17306002	3
Phönix	Ethanol 90% pure / denatured	17306002	1
Caelo	Ethanol 90% pure / denatured	17306205	1
Caelo	Ethanol 90% pure / denatured	17306101	1
Caelo	Ethanol 90% pure / denatured	17306003	1
Tsiripidou	Ethanol 90% pure / denatured	180324	1
Apomix/Noweda	Ethanol 90% pure / denatured		1
Brinn-Holtz	Ethanol 90% pure / denatured	190615	2
Eigenherstellung Hubert...	Ethanol 90% pure / denatured	2014072222DF	1
Caelo	Ethanol 90% pure / denatured	2515A-03042	1
Rosenapotheke	Ethanol 90% pure / denatured	2417E-07094	1
apomix	Ethanol 90% pure / denatured	18A080615	1
Defektur	Ethanol 90% pure / denatured	2115A-03042	1
Kastner (eigen)	Ethanol 90% pure / denatured	3114E-03042	1
Hochhaus-Apotheke	Ethanol 90% pure / denatured	2616Q-03042	1
Apomix / Phoenix	Ethanol 90% pure / denatured	2040030717	1
Eigenherstellung	Ethanol 90% pure / denatured	3416E-07057	1
apomix PKH GmbH Halle	Ethanol 90% pure / denatured	30A310815	1
Eigenherstellung	Ethanol 90% pure / denatured	2614I-02303	1
GEHE	Ethanol 90% pure / denatured	28A110917	1
Defektur	Ethanol 90% pure / denatured	3015M-03042	1
Eigenherstellung	Ethanol 90% pure / denatured	3416E-02303	1
Apotheke	Ethanol 90% pure / denatured	3116Q-03042	1
Eigen/Neue Apotheke	Ethanol 90% pure / denatured	4315I-02303	1
Apomix	Ethanol 90% pure / denatured	37A021214	1
Eigenherstellung	Ethanol 90% pure / denatured	4414A-07057	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apotärze	Ethanol 90% pure / denatured	4414E-03042	1
Engel/Engel	Ethanol 90% pure / denatured	4415A-07014	1
Hochhaus-Apotheke	Ethanol 90% pure / denatured	4514A-03042	1
apomix	Ethanol 90% pure / denatured	45A261015	1
Rondell Apotheke	Ethanol 90% pure / denatured	4615E-03042	1
Apomix	Ethanol 90% pure / denatured	46A071215	1
Rondell Apotheke	Ethanol 90% pure / denatured	4615E-02303	1
Rondell Apotheke	Ethanol 90% pure / denatured	4615E-07209	1
Defektur	Ethanol 90% pure / denatured	4714E-03042	1
Alte Apotheke	Ethanol 90% pure / denatured	4716E-03042	1
Defektur	Ethanol 90% pure / denatured	5015I-03042	1
Apomix / Phoenix	Ethanol 90% pure / denatured	47A071215	1
Defektur	Ethanol 90% pure / denatured	60616/50Hk	1
Rondell Apotheke	Ethanol 90% pure / denatured	4814E-02303	1
Wolf-Apotheke-Defektur	Ethanol 90% pure / denatured	80711/50Hk	1
Wolf Apotheke/Defektur	Ethanol 90% pure / denatured	80615/50Kr	1
Caelo	Ethanol 90% pure / denatured	80933148	1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-160416-1D	1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-030117-1Praxis	1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-021116-3D	1
Tsiripidou (Fortuna Apo...	Ethanol 90% pure / denatured		1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-220716-1P	1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-130216-1D	1
Martinus-Apotheke	Ethanol 90% pure / denatured	E-211014-1D	1
Martinus-Apothke	Ethanol 90% pure / denatured	E-170215-1D	1
	Ethanol 90% pure / denatured	Hergestelllt070514	1
	Ethanol 90% pure / denatured		2
Eigene Herstellung	Ethanol 90% pure / denatured	SP9090	1
Eigene Herstellung	Ethanol 90% pure / denatured	SPV9024	1
Eigenherstellung	Ethanol 90% pure / denatured	4211215	1
Caelo	Ethanol 90% pure / denatured	12275501	1
VDL 09.05.14	Ethanol 90% pure / denatured	13399801	1
Gutenberg-Apotheke	Ethanol 90% pure / denatured	140324001	1
Caelo	Ethanol 90% pure / denatured	17173903	1
Eigenherstellung	Ethanol 90% pure / denatured	1312I-07057	2
Eigenherstellung	Ethanol 90% pure / denatured	514E-07057	2
Eigenherstellung	Ethanol 90% pure / denatured	1613A-07057	1
Eigenherstellung	Ethanol 90% pure / denatured	4313Q-07057	1
Eigenherstellung	Ethanol 90% pure / denatured	4713Q-07057	2
Defektur	Ethanol 90% pure / denatured	26112012D01	1
Apomix	Ethanol 90% pure / denatured	12A150513	1
Defektur	Ethanol 90% pure / denatured	30082012	1
Caelo	Ethanol 90% pure / denatured	814E-02303	1
Nachprüfung	Ethanol 90% pure / denatured	913Q-03042	1
Caelo	Ethanol 90% pure / denatured	915Q-02303	1
Hustadt apotheke	Ethanol 90% pure / denatured	130612/1	1
Post-Apotheke	Ethanol 90% pure / denatured	1414A-02303	1
Caelo	Ethanol 90% pure / denatured	16106801	1
Caelo	Ethanol 90% pure / denatured	1915L-02303	1
Defektur	Ethanol 90% pure / denatured	1816A-03042	1
Defektur	Ethanol 90% pure / denatured	3216N-03042	1
Defektur	Ethanol 90% pure / denatured	4916E-03042	1
Caelo	Ethanol 90% pure / denatured	5214E-02303	1
Wolf-Apotheke / Defektur	Ethanol 90% pure / denatured	70428/50Hk	1
Berkel AHK	Ethanol 90% pure / denatured		1
Caelo	Ethanol 90% pure / denatured	Defektur	1
eigene Herstellung	Ethanol 90% pure / denatured	SP9025	1
Römer-Apotheke	Ethanol 90% pure / denatured	1914E-03042	1
Caelo	Ethanol 90% pure / denatured	14027305	1
Germania Apo	Ethanol 90% pure / denatured	219060	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ethanol 90% pure / denatured	2413Q-03042	1
Römer-Apotheke	Ethanol 90% pure / denatured	36141-03042	1
Rondell-Apotheke	Ethanol 90% pure / denatured	3913Q-03042	1
Römer-Apotheke	Ethanol 90% pure / denatured	4213E-03042	1

- 90 891 spectra from 501 *Apo-Ident* customers from a total of 26 749 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Ethanol 90% pure / denatured* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Ethanol 90% pure / denatured* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	280	0	93 455
Type C	0	456	5	90 782

The substance/substance group *Ethanol 90% pure / denatured* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9765 %)	100.0000 % (> 97.8571 %)
Type C	100.0000 % (> 99.6601 %)	98.9154 % (> 98.2646 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31056	31056	0.00	128.33
31392	31056	11.76	114.40
33102	33102	0.00	104.82

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Ethanol 96% pure / denatured**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31488-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Ethanol 96% pure / denatured; Ethanolum 96% rein / vergällt

Special notes

When selecting the *Ethanol 96% pure / denatured* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Ethanol 96% pure / denatured*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Merck KGaA	Ethanol 96% pure...	K45197671406	31488	40	1503676
Laborchemie A...	Ethanol 96% pure...	21001382	31564	40	20100725

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Ethanol 96% pure / denatured*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 160 spectra from 5 reference samples of the substance/substance group *Ethanol 96% pure / denatured*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Laborchemie Apolda	Ethanol 96% pure / denatured	16000409	30534	40
Merck	Ethanol 96% pure / denatured	K45197671406	31488 [†]	20
Laborchemie Apolda	Ethanol 96% pure / denatured	21001382	31564 [†]	20
Laborchemie Apolda	Ethanol 96% pure / denatured	24001427	33852	40
Laborchemie Apolda	Ethanol 96% pure / denatured	25000305	34123	40

- 93 714 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 664 spectra from 205 *Apo-Ident* customers from 368 batches from the substance/substance group *Ethanol 96% pure / denatured*.
- These include spectra of independent samples from 352 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Alkoholvertrieb Süd GmbH	Ethanol 96% pure / denatured	2/0713	1
Alkoholvertrieb Dr. Ber...	Ethanol 96% pure / denatured	15132014/02122013	1
hetterich	Ethanol 96% pure / denatured	3574702	1
Branntweinmonopol	Ethanol 96% pure / denatured	613Et0H96	1
Applichem/Diagonal	Ethanol 96% pure / denatured	70214PJK05	1
Bundesmonopolverwaltung...	Ethanol 96% pure / denatured		10
Bombastus	Ethanol 96% pure / denatured		1
Caelo	Ethanol 96% pure / denatured	1/12	1
BfB	Ethanol 96% pure / denatured	1/2/3/4	1
Kassner & Sasse	Ethanol 96% pure / denatured	12/11/23	1
Kottkamp/Gehe	Ethanol 96% pure / denatured	121213	1
Bergchemie	Ethanol 96% pure / denatured	13/03/06	2
Bergchemie	Ethanol 96% pure / denatured		2
Berkel AHK	Ethanol 96% pure / denatured	13/07/17-1	1
Berkel AHK	Ethanol 96% pure / denatured	13/11/08	1
Bergchemie	Ethanol 96% pure / denatured	13/11/27	1
Berkel AHK	Ethanol 96% pure / denatured	1410T31040613	1
kottkamp,gehe	Ethanol 96% pure / denatured	130513	1
Alkoholvertrieb Süd GmbH...	Ethanol 96% pure / denatured	1313081029/26082013	1
bergchemie	Ethanol 96% pure / denatured	14-06-04	1
Berkel AHK	Ethanol 96% pure / denatured	1410T31191113	1
Bergchemie	Ethanol 96% pure / denatured	14-08-12	2
Berkel AHK	Ethanol 96% pure / denatured	1410T311060513	1
Tereos 20131470D	Ethanol 96% pure / denatured	1410T31040613	1
Miller	Ethanol 96% pure / denatured	1413040904/16042013	1
Noweda	Ethanol 96% pure / denatured	19000663	2
Apolda GmbH	Ethanol 96% pure / denatured	19001936	1
Laborchemi Apolda	Ethanol 96% pure / denatured	19002423	1
BMV Branntwein	Ethanol 96% pure / denatured		12
Gehe 18.04.2013	Ethanol 96% pure / denatured	20000065	1
Laborchemie Apolda	Ethanol 96% pure / denatured	20000065	2
Laborchemie/Gehe	Ethanol 96% pure / denatured	20000071	1
Laborchemie Apolda	Ethanol 96% pure / denatured	20000071	1
Laborchemie/Gehe	Ethanol 96% pure / denatured	20000305	3
Apolda GmbH	Ethanol 96% pure / denatured	20000305	1
Heyl	Ethanol 96% pure / denatured	20000947	1
Noweda / Apolda	Ethanol 96% pure / denatured	20000974	2
Gehe/Laborchemie Apolda	Ethanol 96% pure / denatured	20000973	2
Apolda / Gehe	Ethanol 96% pure / denatured	20000974	1
Laborchemie Apolda/Noweda	Ethanol 96% pure / denatured	20000974	1
Apolda, Noweda	Ethanol 96% pure / denatured	20001212	1
Laborchemie APOLDA/ Now...	Ethanol 96% pure / denatured	20001574	1
Noweda/Labochemie	Ethanol 96% pure / denatured	20001566	2
Sanacorp	Ethanol 96% pure / denatured	20001574	1
Laborchemie Apolda	Ethanol 96% pure / denatured	20001886	2
Apolda/VDL	Ethanol 96% pure / denatured	20001830	1
Spangro/Laborchemie Apo...	Ethanol 96% pure / denatured	20001583	1
Heyl / AHD	Ethanol 96% pure / denatured	20001830	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Laborchemie Apolda	Ethanol 96% pure / denatured	20002191	1
Gehe am 18.02.2014	Ethanol 96% pure / denatured	20002191	1
Laborchemie Apolda Heyl...	Ethanol 96% pure / denatured	20001887	1
Apolda / Sanacorp	Ethanol 96% pure / denatured	20002191	1
W. Hoffmann/Noweda	Ethanol 96% pure / denatured	20120106	2
Schmittmann GmbH	Ethanol 96% pure / denatured	20130105	1
Schmittmann GmbH	Ethanol 96% pure / denatured	20130105	1
Bundesmonopolverw.	Ethanol 96% pure / denatured	20130903	2
Laborchemie APOLDA	Ethanol 96% pure / denatured	21002269	1
Laborchemie Apolda	Ethanol 96% pure / denatured	22001720	1
Caelo	Ethanol 96% pure / denatured	2310A0341	1
Werner Hofmann	Ethanol 96% pure / denatured	30091302	1
Apolda/Phönix	Ethanol 96% pure / denatured	3613M-02304	1
Berkel AHK	Ethanol 96% pure / denatured	3613Q2304	1
Kost/kost/sophien	Ethanol 96% pure / denatured	4065280513	1
kost	Ethanol 96% pure / denatured	410031212052	1
Kost/Krieger/04.07.13	Ethanol 96% pure / denatured	410040713061	2
Kost/13.12.13	Ethanol 96% pure / denatured	410051213065	4
kost	Ethanol 96% pure / denatured	410031212057	1
Sophien-Arzneimittel/dir...	Ethanol 96% pure / denatured	410081013063	2
Kost/kost/sophien	Ethanol 96% pure / denatured	410051213065	2
Kost/Krieger/14.10.13	Ethanol 96% pure / denatured	410090813062	3
Kost/Krieger/05.12.13	Ethanol 96% pure / denatured	410111113064	3
Kost/Krieger/12.08.13	Ethanol 96% pure / denatured	410130313059	2
Kost/Krieger/26.06.13	Ethanol 96% pure / denatured	410140513060	1
Kost/Krieger/26.06.13	Ethanol 96% pure / denatured	410150513060	1
Kost/Krieger/21.02.14	Ethanol 96% pure / denatured	410220114066	3
Sonnenberg/Zentrum	Ethanol 96% pure / denatured	4711E-03041	1
Brüggenmann/Alkohol	Ethanol 96% pure / denatured	7814060362	1
BMV Branntwein	Ethanol 96% pure / denatured	Ch19	1
Werner Hoffmann	Ethanol 96% pure / denatured	Ch20130033	1
Hoffmann	Ethanol 96% pure / denatured	E-200412-1	4
Merck	Ethanol 96% pure / denatured	K43729871232	3
Brüggenmann/Alkohol	Ethanol 96% pure / denatured	7814020283	1
Monopolstelle Jenne	Ethanol 96% pure / denatured	J150114	1
Hoffmann/ Noweda	Ethanol 96% pure / denatured	E-231213-1	3
Merck	Ethanol 96% pure / denatured	K44538471324	1
BFB	Ethanol 96% pure / denatured	keinevorhanden	2
Gehe	Ethanol 96% pure / denatured	LieferungMai2013	1
Taoasis/K	Ethanol 96% pure / denatured	VE30892-112552	1
BfB-Betrieb München	Ethanol 96% pure / denatured		1
bergchemie	Ethanol 96% pure / denatured	UN1170	1
PVL	Ethanol 96% pure / denatured	1204131	1
Bergchemie	Ethanol 96% pure / denatured	10/01/10	1
Werner Hofmann	Ethanol 96% pure / denatured	8021304	1
Bergchemie	Ethanol 96% pure / denatured	12/06/14	1
Berkel AHG	Ethanol 96% pure / denatured	12/11/23	1
Noweda	Ethanol 96% pure / denatured	1202124	2
Caelo	Ethanol 96% pure / denatured	12345	1
Klenk	Ethanol 96% pure / denatured	123456	1
Phönix/berkel	Ethanol 96% pure / denatured	130208	1
Herbasin	Ethanol 96% pure / denatured	13/01/10	1
Berkel AHK	Ethanol 96% pure / denatured	1410T31141212	1
Berkel AHK	Ethanol 96% pure / denatured	1410T31080113	1
Berkel AHK	Ethanol 96% pure / denatured	170412/292	1
Euro OTC	Ethanol 96% pure / denatured	18001042	1
Laborch. Apolda	Ethanol 96% pure / denatured	18001881	1
Anzag	Ethanol 96% pure / denatured	19000121	1
Laborchemie Apolda GmbH...	Ethanol 96% pure / denatured	19000121	1
Gehe Laborchemie Apolda...	Ethanol 96% pure / denatured	19000689	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Anzag 24.05.2012	Ethanol 96% pure / denatured	19000700	1
Laborchemie APOLDA GmgH	Ethanol 96% pure / denatured	19000689	1
Laborchemie Apolda	Ethanol 96% pure / denatured	19001519	2
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	19001937	1
Noweda	Ethanol 96% pure / denatured	19002422	2
Apolda / Gehe	Ethanol 96% pure / denatured	20000065	1
Schmittmann/Hoffmann	Ethanol 96% pure / denatured	20020112-1	1
Bioenergie Icking	Ethanol 96% pure / denatured	201102170016	1
Werner Hofmann	Ethanol 96% pure / denatured	20150278	1
Hoffmann	Ethanol 96% pure / denatured	28081201	1
Berkel AHK	Ethanol 96% pure / denatured	2906211085	1
Applichem	Ethanol 96% pure / denatured	220113PML01	1
Caelo	Ethanol 96% pure / denatured		3
Fagron	Ethanol 96% pure / denatured		1
Werner Hoffmann	Ethanol 96% pure / denatured		1
Caelo	Ethanol 96% pure / denatured	5678	1
Gehe	Ethanol 96% pure / denatured	DE00001237810	1
sanacorp	Ethanol 96% pure / denatured	K40787471	1
VDL	Ethanol 96% pure / denatured	K43005871151	2
Caelo	Ethanol 96% pure / denatured	K43619871227	1
Berkel AHK	Ethanol 96% pure / denatured	VE28587	1
Noweda	Ethanol 96% pure / denatured	19001935	1
Hersteller: Laborchemie...	Ethanol 96% pure / denatured	172N0346	1
Berkel	Ethanol 96% pure / denatured	1410T31260814	1
Apolda Laborchemie	Ethanol 96% pure / denatured	19002599	1
DHL	Ethanol 96% pure / denatured	11062013	1
Werner Hofmann	Ethanol 96% pure / denatured	111111	1
Alkoholvertrieb Süd GmbH	Ethanol 96% pure / denatured	1313081029/26082013	1
Phoenix	Ethanol 96% pure / denatured	200000973	1
Laborchemie Apolda	Ethanol 96% pure / denatured	20001199	1
Laborchemie / Noweda	Ethanol 96% pure / denatured	20001830	1
Apolda Laborchemie	Ethanol 96% pure / denatured	20001887	1
Laborchemie Apolda	Ethanol 96% pure / denatured	20001885	1
Laborchemie	Ethanol 96% pure / denatured	21000152	1
Apolda GmbH	Ethanol 96% pure / denatured	22001717	1
Brand GmbH Alkoholvertr...	Ethanol 96% pure / denatured	24052013/60	1
spinnrad	Ethanol 96% pure / denatured	815310072	1
Alkoholvertrieb Süd GmbH	Ethanol 96% pure / denatured	keine	1
Kammer Kirschner	Ethanol 96% pure / denatured		3
Alkoholvertrieb Süd GmbH...	Ethanol 96% pure / denatured	6282017/24102016	2
Alkoholvertrieb Sued	Ethanol 96% pure / denatured	87872017/29052017	1
Agralko	Ethanol 96% pure / denatured	1062015	1
Phönix	Ethanol 96% pure / denatured	18122017/25112016	1
Schmittmann GmbH	Ethanol 96% pure / denatured	2041503	1
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	9012016/05112015	1
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	18122017/25112016	1
Allianz Healthcare	Ethanol 96% pure / denatured	2394894	2
Tereos/Berkel AHK	Ethanol 96% pure / denatured	31115/705	1
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	3297636/11062014	1
kottkamp/Anzag	Ethanol 96% pure / denatured	61015/505/017	1
Alkoholvertrieb Süd/Gehe	Ethanol 96% pure / denatured	3297636/11062014	1
Senge& Butz GmbH WE: 0...	Ethanol 96% pure / denatured	41209	1
Alkoholvertrieb Süd Dr. ...	Ethanol 96% pure / denatured	39552016/28012016	1
Werner Hoffmann	Ethanol 96% pure / denatured	7081503	1
Alkoholvertrieb süd	Ethanol 96% pure / denatured	53472018/0032018	1
Alkoholvertrieb süd	Ethanol 96% pure / denatured	71192018/23042018	1
Apolda GmbH	Ethanol 96% pure / denatured	100050	1
Berkel AHK	Ethanol 96% pure / denatured	40517/843	1
Berkel AHK	Ethanol 96% pure / denatured	101017/884	1
Bundes...	Ethanol 96% pure / denatured	8042014	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	80452016/01062016	1
Berkel AHK	Ethanol 96% pure / denatured	111016/786	2
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	116702016/05092016	1
Berkel AHK	Ethanol 96% pure / denatured	120516/746	1
Alkohol Vertrieb Süd	Ethanol 96% pure / denatured	121952017/21082017	1
Taoasis	Ethanol 96% pure / denatured	122428	1
Laborchemie Apolda/Jenn...	Ethanol 96% pure / denatured	12345	4
Heyl / Sanacorp	Ethanol 96% pure / denatured		1
Hoffmann	Ethanol 96% pure / denatured	123456	1
Kottkamp/Anzag	Ethanol 96% pure / denatured	130116/505/020	1
Berkel AHK	Ethanol 96% pure / denatured	14/01/30	1
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	138512107/28092017	1
Bergchemie	Ethanol 96% pure / denatured	14/02/27	1
Berkel AHK	Ethanol 96% pure / denatured	14/10/16	1
Berkel AHK	Ethanol 96% pure / denatured	14/09/17	1
Berkel AHK	Ethanol 96% pure / denatured	14/11/13	1
/Phönix	Ethanol 96% pure / denatured	140604	1
Bergchemie	Ethanol 96% pure / denatured	14-06-04	1
Bundesm.Brandwein	Ethanol 96% pure / denatured	140710	1
Bergchemie	Ethanol 96% pure / denatured	14-10-16	1
Berkel AHK	Ethanol 96% pure / denatured	1410T31050914	1
Berkel	Ethanol 96% pure / denatured	1410T31051115	1
Berkel	Ethanol 96% pure / denatured	1410T31141015	1
Fagron	Ethanol 96% pure / denatured	1410T31200516	1
Berkel AHK	Ethanol 96% pure / denatured	1410T31280116	1
BerkelAHK	Ethanol 96% pure / denatured	1410T31271115	1
bergchemie	Ethanol 96% pure / denatured	14-12-10	1
Berkel	Ethanol 96% pure / denatured	1410T31310715	1
Bergchemie	Ethanol 96% pure / denatured	15/06/10	1
Berkel AHK	Ethanol 96% pure / denatured	15/03/15	1
Berkel	Ethanol 96% pure / denatured	15/06/10	1
Berkel AHK	Ethanol 96% pure / denatured	15/07/23	2
Bergchemie	Ethanol 96% pure / denatured	15-05-07	1
Bergchemie	Ethanol 96% pure / denatured	15-10-09	1
Berkel AHK	Ethanol 96% pure / denatured	151125	1
Kottkamp Alkoholvertrie...	Ethanol 96% pure / denatured	1512155050519	1
Bergchemie	Ethanol 96% pure / denatured	15-5-07	1
Bergchemie	Ethanol 96% pure / denatured	16/01/06	2
Bergchemie	Ethanol 96% pure / denatured	16/08/03	1
Bergchemie	Ethanol 96% pure / denatured	16/02/02	1
BerkelAHK	Ethanol 96% pure / denatured	16/11/23	1
Kottkamp / über GEHE	Ethanol 96% pure / denatured	160714/505/004	1
Brüggemann	Ethanol 96% pure / denatured	1514041040	1
Hoffmann	Ethanol 96% pure / denatured	16061405	1
VWR	Ethanol 96% pure / denatured	16E034127	1
Laborchemie Apolda	Ethanol 96% pure / denatured	17041401	1
Bergchemie/Bergchemie	Ethanol 96% pure / denatured	17-05-19	2
Berkel AHK	Ethanol 96% pure / denatured	180327	1
Fagron	Ethanol 96% pure / denatured	19001836	1
Caelo	Ethanol 96% pure / denatured	20001583	1
Laborchemie/Noweda	Ethanol 96% pure / denatured	19001836	1
Gehe-Apolda	Ethanol 96% pure / denatured	2000/1895	1
Gehe-Apolda	Ethanol 96% pure / denatured	20001885	2
Apolda/Gehe	Ethanol 96% pure / denatured	20002193	1
Kassner & Sasse	Ethanol 96% pure / denatured	200-578-6	1
Eal Euro-Alkohol	Ethanol 96% pure / denatured	20140100	1
Schmittmann GmbH	Ethanol 96% pure / denatured	20140100	1
Schmittmann	Ethanol 96% pure / denatured	20150108	1
Hoffmann/Schmittmann	Ethanol 96% pure / denatured	20150108	1
Agralko	Ethanol 96% pure / denatured	2015010901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro Alkohol GmbH	Ethanol 96% pure / denatured	20150108	1
Caelo	Ethanol 96% pure / denatured	20150908	1
eal/euro-alkohol	Ethanol 96% pure / denatured	20150108	1
Hofmann/Schmittmann	Ethanol 96% pure / denatured	20160007	1
Werner Hofmann	Ethanol 96% pure / denatured	20170160	1
Laborchemie Apolda	Ethanol 96% pure / denatured	21000152	2
Laborchemie Apolda/Phönix	Ethanol 96% pure / denatured	21000152	1
Noweda	Ethanol 96% pure / denatured	21000197	1
Apolda/ Sanacorp	Ethanol 96% pure / denatured	21000154	2
LaborchemieApolda/Sanac...	Ethanol 96% pure / denatured	21001053	1
Noweda / Aploda	Ethanol 96% pure / denatured	21000197	1
Apolda/Phoenix	Ethanol 96% pure / denatured	21000197	1
All.Heal./laborchemie A...	Ethanol 96% pure / denatured	21001053	1
Heyl / Sanacorp	Ethanol 96% pure / denatured	21001373	1
Apolda GmbH	Ethanol 96% pure / denatured	21001379	2
Laborchemie Apolda	Ethanol 96% pure / denatured	21001054	1
alliance Health. 23.10.14	Ethanol 96% pure / denatured	21001379	1
Apolda/ Sanacorp	Ethanol 96% pure / denatured	21001560	1
Apolda/ Sanacorp	Ethanol 96% pure / denatured	21001380	1
Apolda GmbH	Ethanol 96% pure / denatured	21001382	1
Heyl / Sanacorp	Ethanol 96% pure / denatured	21001379	1
All.Heal./laborchemie A...	Ethanol 96% pure / denatured	21001561	1
Phoenix	Ethanol 96% pure / denatured	21001561	1
Laborchemie Apolda über...	Ethanol 96% pure / denatured	21001964	1
Caelo	Ethanol 96% pure / denatured	21001964	1
Apolda gmbH/N	Ethanol 96% pure / denatured	21002064	1
Laborchemie Apolda	Ethanol 96% pure / denatured	21002075	3
AHD 03.12.14 EK 85,21	Ethanol 96% pure / denatured	21002075	1
Heyl(Laborchemie Apolda...	Ethanol 96% pure / denatured	21002064	1
Laborchemie Apolda/Phoe...	Ethanol 96% pure / denatured	21001736	1
Phönix 06.2.15	Ethanol 96% pure / denatured	21002075	1
Gehe/Dr. Wollff	Ethanol 96% pure / denatured	21002237	2
Noweda/laborchemie Apolda	Ethanol 96% pure / denatured	21002237	1
Laborch. Apolda/ Noweda	Ethanol 96% pure / denatured	21002237	1
Caelo	Ethanol 96% pure / denatured	21002244	1
Laborchemie Apolda	Ethanol 96% pure / denatured	21002244	2
Apolda/ Sanacorp	Ethanol 96% pure / denatured	21002239	1
Sanacorp	Ethanol 96% pure / denatured	21002268	1
Apolda	Ethanol 96% pure / denatured	21002268	1
Laborchemie Apolda/Phönix	Ethanol 96% pure / denatured	21002268	1
Apolda/Phönix	Ethanol 96% pure / denatured	21002269	1
Kottkamp Alkoholvertrie...	Ethanol 96% pure / denatured	211114/505/008	1
Caelo	Ethanol 96% pure / denatured	21002269	1
Laborchemie Apold GmbH	Ethanol 96% pure / denatured	21002269	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	21002268	1
Caelo	Ethanol 96% pure / denatured	220000086	1
Laborchemie Apolda/Noweda	Ethanol 96% pure / denatured	22000086	1
Phoenix Goettin	Ethanol 96% pure / denatured	21002269	1
Apolda/ Sanacorp	Ethanol 96% pure / denatured	22000090	1
Apolda/Sanacorp	Ethanol 96% pure / denatured	22000090	1
Caelo	Ethanol 96% pure / denatured	22000096	1
Laborchemie/Noweda	Ethanol 96% pure / denatured	22000161	1
Fagron	Ethanol 96% pure / denatured	22000096	1
Laborchemie Apolda/Alli...	Ethanol 96% pure / denatured	22000392	1
Heyl	Ethanol 96% pure / denatured	22000392	1
Heyl/Noweda	Ethanol 96% pure / denatured	22000326	1
Laborchemie Apolda/Gehe	Ethanol 96% pure / denatured	22000394	2
Caelo	Ethanol 96% pure / denatured	22000392	1
Laborchemie Apolda/Phönix	Ethanol 96% pure / denatured	22000392	1
Heyl	Ethanol 96% pure / denatured	22000394	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Laborchemie/Gehe	Ethanol 96% pure / denatured	22000603	3
Laborchemie Apolda	Ethanol 96% pure / denatured	22000394	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	22000394	1
Laborchemie Apolda GmbH.	Ethanol 96% pure / denatured	22000603	1
Apolda/Phönix	Ethanol 96% pure / denatured	22000616	1
Laborchemie Apolda GmbH...	Ethanol 96% pure / denatured	22000616	1
Laborchemie Apolda	Ethanol 96% pure / denatured	22000616	1
Laborchemie/Sanacorp	Ethanol 96% pure / denatured	22000616	1
Noweda 04.08.15	Ethanol 96% pure / denatured	22000616	1
Caelo	Ethanol 96% pure / denatured	22000641	1
Apolda/Phoenix	Ethanol 96% pure / denatured	22000616	1
Laborchemie Apolda	Ethanol 96% pure / denatured	22000999	1
All. Healthc. / Laborph...	Ethanol 96% pure / denatured	22000616	1
Caelo	Ethanol 96% pure / denatured	22000999	1
Laborchemie Apolda/Sana...	Ethanol 96% pure / denatured	22000999	1
Apolda gmbH/N	Ethanol 96% pure / denatured	22000999	1
Apolda/Phönix	Ethanol 96% pure / denatured	22001291	1
Laborchemie Apolda/Alli...	Ethanol 96% pure / denatured	22001291	1
Laborchemie Apolda/Phönix	Ethanol 96% pure / denatured	22001291	1
heyl	Ethanol 96% pure / denatured	22001291	1
Laborchemie Apolda	Ethanol 96% pure / denatured	22001292	1
Laborchemie Apolda	Ethanol 96% pure / denatured	22001385	1
AllianceHealthcare	Ethanol 96% pure / denatured	22001291	1
Laborchemie / Sanacorp	Ethanol 96% pure / denatured	22001717	1
Laborch. Apolda/ Noweda	Ethanol 96% pure / denatured	22001291	1
Heyl/Gehe	Ethanol 96% pure / denatured	22001717	1
Laborch. Apolda/ Noweda	Ethanol 96% pure / denatured	22001718	1
Apolda / Anzag	Ethanol 96% pure / denatured	22001717	1
Apolda	Ethanol 96% pure / denatured	22001718	1
Hedinger	Ethanol 96% pure / denatured	22001720	1
Caelo	Ethanol 96% pure / denatured	22001718	1
Bionorica	Ethanol 96% pure / denatured	22001720	1
Laborchemie Apold GmbH	Ethanol 96% pure / denatured	22001948	1
Apolda GmbH	Ethanol 96% pure / denatured	22001948	1
Laborchemie Apolda/Phoe...	Ethanol 96% pure / denatured	22001948	2
Heyl/Gehe	Ethanol 96% pure / denatured	22001948	1
Apolda / Sanacorp	Ethanol 96% pure / denatured	22001985	1
Caelo	Ethanol 96% pure / denatured	22001985	1
Apolda	Ethanol 96% pure / denatured	22002230	1
Laborchemie Apolda über...	Ethanol 96% pure / denatured	22002230	1
Laborchemie Apolda/ San...	Ethanol 96% pure / denatured	22002231	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	22001985	1
Laborchemie Apoldo GmbH	Ethanol 96% pure / denatured	22002231	1
Caelo	Ethanol 96% pure / denatured	23000175	1
Gehe	Ethanol 96% pure / denatured	23000175	1
Heyl	Ethanol 96% pure / denatured	23000188	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	23000188	2
Laborchemie Apolda	Ethanol 96% pure / denatured	23000188	1
Heyl	Ethanol 96% pure / denatured	23000189	1
Caelo	Ethanol 96% pure / denatured	23000480	1
Apolda/Phönix	Ethanol 96% pure / denatured	23000749	2
Laborchemie Apolda/ Now...	Ethanol 96% pure / denatured	23000750	6
Caelo	Ethanol 96% pure / denatured	23000749	1
Heyl/Noweda	Ethanol 96% pure / denatured	23000749	1
LABORCHEMIE/Gehe	Ethanol 96% pure / denatured	23000750	1
Sanacop	Ethanol 96% pure / denatured	23000780	1
Laborchemie APOLDA GmbH	Ethanol 96% pure / denatured	23000997	1
Laborch. Apolda/ Noweda	Ethanol 96% pure / denatured	23000997	1
Laborchemie Apolda über...	Ethanol 96% pure / denatured	23000997	1
AllianceHealthcare	Ethanol 96% pure / denatured	23000997	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Laborchemie Apolda	Ethanol 96% pure / denatured	23000997	1
Laborchemie Apolda/ Now...	Ethanol 96% pure / denatured	23001055	4
Laborchemie Apolda/Noweda	Ethanol 96% pure / denatured	23000997	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	23001368	1
Laborchemie Apolda GmbH...	Ethanol 96% pure / denatured	23001357	1
Laborchemie Apolda/Gehe	Ethanol 96% pure / denatured	23001357	1
Apolda / Gehe	Ethanol 96% pure / denatured	23001357	1
Caelo	Ethanol 96% pure / denatured	23001357	1
Caelo	Ethanol 96% pure / denatured	23001368	2
Laborchemie Apolda/ All...	Ethanol 96% pure / denatured	23001368	1
Apolda	Ethanol 96% pure / denatured	23001368	1
Laborchemie Apolda/Phönix	Ethanol 96% pure / denatured	23001368	1
Heyl / Noweda	Ethanol 96% pure / denatured	23001368	1
Heyl/Phönix	Ethanol 96% pure / denatured	23001368	1
laborchemie apolda gmbh	Ethanol 96% pure / denatured	23001368	1
Laborchemie Apolda GmbH...	Ethanol 96% pure / denatured	23001368	1
Apolda/Noweda	Ethanol 96% pure / denatured	23001368	1
Laborchemie Apolda/ Now...	Ethanol 96% pure / denatured	23001785	4
Gehe	Ethanol 96% pure / denatured	23001968	1
Heyl Laborchemie Apolda...	Ethanol 96% pure / denatured	23002116	1
Laborchemie Apold GmbH	Ethanol 96% pure / denatured	23001968	1
Heyl/Noweda	Ethanol 96% pure / denatured	23002116	2
LABORCHEMIE APOLDA GmbH...	Ethanol 96% pure / denatured	23001968	1
Hubertus Apotheke	Ethanol 96% pure / denatured	23001968	1
Fagron	Ethanol 96% pure / denatured	23002116	1
Kottkamp Alkoholvertr.	Ethanol 96% pure / denatured	230414	1
Apolda/Noweda	Ethanol 96% pure / denatured	24000163	1
Laborchemie Apolda	Ethanol 96% pure / denatured	24000161	1
Laborchemie Apolda Heyl...	Ethanol 96% pure / denatured	24000186	1
LABORCHEMIE/Gehe	Ethanol 96% pure / denatured	24000163	1
Laborchemie Apolda	Ethanol 96% pure / denatured	24000188	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	24000186	1
Phönix/Heyl	Ethanol 96% pure / denatured	24000863	1
Phönix / Heyl	Ethanol 96% pure / denatured	24000856	1
LABORCHEMIE/Gehe	Ethanol 96% pure / denatured	24000856	1
Sanacorp	Ethanol 96% pure / denatured	24001427	1
Laborchemie Apolda	Ethanol 96% pure / denatured	24001427	1
Apolda/Sanacorp	Ethanol 96% pure / denatured	24002163	1
Apolda GmbH	Ethanol 96% pure / denatured	24002163	1
Euro OTC	Ethanol 96% pure / denatured	27042017	1
Heyl	Ethanol 96% pure / denatured	25000305	1
anzag	Ethanol 96% pure / denatured	25439	1
Berkel AHK	Ethanol 96% pure / denatured	260411/145	2
Berkel AHK	Ethanol 96% pure / denatured	310317/835	1
Hessische Spezialitäten...	Ethanol 96% pure / denatured	4101081501	1
Kost/sophien	Ethanol 96% pure / denatured	410010414067	1
Kwst/Walther CMP	Ethanol 96% pure / denatured	290814	1
Rondell Apotheke	Ethanol 96% pure / denatured	410040216085	1
Fa. Karl Kost Koblenz	Ethanol 96% pure / denatured	410041115082	1
Fa. Karl Kost Koblenz	Ethanol 96% pure / denatured	410040216085	1
Miller Alkoholvertrieb ...	Ethanol 96% pure / denatured	410-103-150508/26052	1
Karl-Josef Kost Alkohol...	Ethanol 96% pure / denatured	41004115082	1
Alkoholvertrieb Süd/DPD	Ethanol 96% pure / denatured	41010315061020072015	1
Karl-Josef KG	Ethanol 96% pure / denatured	410110117094	1
Karl-Josef Kost	Ethanol 96% pure / denatured	410181217102	1
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	410-104-1405/2909201	1
Karl-Josef Kost Alkohol...	Ethanol 96% pure / denatured	410190516087	1
Alkoholvertrieb Süd	Ethanol 96% pure / denatured	5643/VALU38679/23	1
Merck	Ethanol 96% pure / denatured	45056971	1
SAV LP	Ethanol 96% pure / denatured	601826	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Kottkamp	Ethanol 96% pure / denatured	71026022	1
Brüggemann Alcohol	Ethanol 96% pure / denatured	7814070382	1
Brüggemann Alcohol	Ethanol 96% pure / denatured	7814120470	1
Brüggemann	Ethanol 96% pure / denatured	7815030598	1
Brüggemann	Ethanol 96% pure / denatured	7814110456	1
Brüggemann Alkohol	Ethanol 96% pure / denatured	7815040634	1
Brüggemann	Ethanol 96% pure / denatured	7815040634	1
Brüggemann Alcohol	Ethanol 96% pure / denatured	7815060735	2
Brüggemann	Ethanol 96% pure / denatured	7815080789	3
Brüggemann	Ethanol 96% pure / denatured	7815110957	1
Brüggemann	Ethanol 96% pure / denatured	7816041610	1
Caelo	Ethanol 96% pure / denatured	90309001/PR	1
Bergchemie	Ethanol 96% pure / denatured	aaaaaa	1
Kassner & Sasse/Apothek...	Ethanol 96% pure / denatured	A2013000183	2
AgrAlko	Ethanol 96% pure / denatured	A2014000328	2
Agraralkohol AG München	Ethanol 96% pure / denatured		1
Alfred Quadflieg, Gelse...	Ethanol 96% pure / denatured	K16/08/03	1
Merck	Ethanol 96% pure / denatured	K45475971	1
Hoffmann/ Martinus-Apot...	Ethanol 96% pure / denatured	E-060815-1	1
Merck	Ethanol 96% pure / denatured	Index-No603-002-00-5	1
Merck	Ethanol 96% pure / denatured	K45475971416	3
Caelo	Ethanol 96% pure / denatured	K46109471	1
Bioenergi Icking GmbH	Ethanol 96% pure / denatured	N781601635	8
Merck	Ethanol 96% pure / denatured	K46405471507	2
Brüggemann	Ethanol 96% pure / denatured	N781601257	1
Brüggemann Alcohol	Ethanol 96% pure / denatured	N781601503	1
Bioenergi Icking GmbH	Ethanol 96% pure / denatured	N781062528	1
Bioenergi Icking GmbH	Ethanol 96% pure / denatured	N781602528	3
BrüggemannAlcohol Witte...	Ethanol 96% pure / denatured	N781703809	5
Brüggemann	Ethanol 96% pure / denatured	N781704589	1
Bioenergi Icking GmbH	Ethanol 96% pure / denatured	N781704668	4
BrüggemannAlcohol Witte...	Ethanol 96% pure / denatured	N781704731	4
Brüggemann	Ethanol 96% pure / denatured	n781704927	1
Brüggemann	Ethanol 96% pure / denatured	V781704599(1)	1
Brüggemann	Ethanol 96% pure / denatured	V781704599(3)	1
Brüggemann	Ethanol 96% pure / denatured	V781704599(4)	1
Brüggemann	Ethanol 96% pure / denatured	V781704599(2)	1
Taoasis / G	Ethanol 96% pure / denatured	VE31397-114015	1
Berkel AHK Alkoholhandel	Ethanol 96% pure / denatured	VE32323	1
Berkel AHK Alkoholhandel	Ethanol 96% pure / denatured	VE32507	1
Berkel AHK	Ethanol 96% pure / denatured	VE37799	1
Taoasis	Ethanol 96% pure / denatured	VE36425-122428	1
Gehe Kassel	Ethanol 96% pure / denatured		1
Gutenberg-Apotheke	Ethanol 96% pure / denatured	904214	1
Kottkamp/Phönix	Ethanol 96% pure / denatured	151215/505/019	1
Kottkamp Alkoholvertrie...	Ethanol 96% pure / denatured	160714/505/004	1
Laborchemie Apolda	Ethanol 96% pure / denatured	21000197	1
Laborchemie APOLDA GmbH	Ethanol 96% pure / denatured	21001561	1
Apolda, direkt best.	Ethanol 96% pure / denatured	21001560	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	21001736	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	21002269	1
Laborchemie Apolda ,Gehe	Ethanol 96% pure / denatured	21002075	1
Apolda/Phönix	Ethanol 96% pure / denatured	22000096	1
Laborchemie APOLDA GmbH	Ethanol 96% pure / denatured	22000616	1
Laborchemie Apolda ,Gehe	Ethanol 96% pure / denatured	22000641	1
Laborchemie APOLDA GmbH	Ethanol 96% pure / denatured	22001291	1
Laborchemie Apolda GmbH	Ethanol 96% pure / denatured	22002231	1
Apolda/Spangro	Ethanol 96% pure / denatured	22002231	1
.....	Ethanol 96% pure / denatured	22001720	1
Anzag/Heyl	Ethanol 96% pure / denatured	25000305	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alkoholvertrieb Dr. Ber...	Ethanol 96% pure / denatured	41010414052992014	1
Caelo	Ethanol 96% pure / denatured	19002599	1
Apolda, Phönix	Ethanol 96% pure / denatured	20000307	1
Apolda GmbH/N/PZN:10137...	Ethanol 96% pure / denatured	20000975	1
A.H.D.	Ethanol 96% pure / denatured	20000305	1
Apolda, Gehe	Ethanol 96% pure / denatured	20001887	1
Sanacorp	Ethanol 96% pure / denatured	20001830	1
Bundesmonopolstelle für...	Ethanol 96% pure / denatured	xxxx	2
Apolda, Gehe	Ethanol 96% pure / denatured	20002192	1
	Ethanol 96% pure / denatured		1
Bundesmonopolverwaltung...	Ethanol 96% pure / denatured	2032012GI	2
Caelo	Ethanol 96% pure / denatured	7813010053	2
Fiebig	Ethanol 96% pure / denatured	12020161/1	1
Kottkamp Alkoholvertrieb	Ethanol 96% pure / denatured	180213	1
Brüggemann	Ethanol 96% pure / denatured	7813050119	1
Brüggemann	Ethanol 96% pure / denatured	7813110228	1
Bergchemie	Ethanol 96% pure / denatured	5/0810	1
Berkel AHK	Ethanol 96% pure / denatured	80318/931	1
Berkel AHK	Ethanol 96% pure / denatured	91214/632	1
Caelo	Ethanol 96% pure / denatured	915A-03041	1
Berkel AHK	Ethanol 96% pure / denatured	1137111/08	1
Post-Apotheke	Ethanol 96% pure / denatured	111a-02304	1
Brüggemann	Ethanol 96% pure / denatured	1514020991	1
Finteler Apotheke	Ethanol 96% pure / denatured	170407	1
BERKEL AHK GmbH	Ethanol 96% pure / denatured	170712/334	4
Laborchemie Apolda/Jenn...	Ethanol 96% pure / denatured	20002191	1
Berkel AHK /TEREOS	Ethanol 96% pure / denatured	190515/669	1
Lamotte/PX	Ethanol 96% pure / denatured	200201121	1
Hofmann,Branntweinmonopol	Ethanol 96% pure / denatured	20140162	1
Schmittmann GmbH	Ethanol 96% pure / denatured	20160426	1
eal	Ethanol 96% pure / denatured	20150108	1
Berkel AHK	Ethanol 96% pure / denatured	200214/567	1
Laborchemie Apolda	Ethanol 96% pure / denatured	21001380	1
Hofmann	Ethanol 96% pure / denatured	20150278	1
Berkel AHK	Ethanol 96% pure / denatured	310317/834	1
Laborchemi Apolda	Ethanol 96% pure / denatured	22000395	1
Berkel AHK	Ethanol 96% pure / denatured	210714/598	1
Laborchemie Apolda GMBH...	Ethanol 96% pure / denatured	21002244	1
Anzag/Hedinger	Ethanol 96% pure / denatured	23000749	1
Karl-Josef Kost	Ethanol 96% pure / denatured	410040216085	1
AgrAlko	Ethanol 96% pure / denatured	410-1405	1
Sophien-AM / K.J.Kost	Ethanol 96% pure / denatured	410041115082	1
Berkel AHK	Ethanol 96% pure / denatured	60919/1La	1
Merck	Ethanol 96% pure / denatured		1
Spinnrad	Ethanol 96% pure / denatured	713310072	1
Sanacorp	Ethanol 96% pure / denatured		1
Brüggemann	Ethanol 96% pure / denatured	7814060362	1
Hofmann,Branntweinmonopol	Ethanol 96% pure / denatured		2
BfB	Ethanol 96% pure / denatured	10065927	2
Berkel AHK	Ethanol 96% pure / denatured	13/09/05	2
Berkel AHK	Ethanol 96% pure / denatured	1410T31181213	1
Berkel	Ethanol 96% pure / denatured	1410T31191113	1
Caelo	Ethanol 96% pure / denatured	20001887	1
Werner Hofmann	Ethanol 96% pure / denatured	20120106	2
Werner Hoffmann	Ethanol 96% pure / denatured	20130033	1
Caelo	Ethanol 96% pure / denatured	AKD7	2
eal Euro-Alkohol	Ethanol 96% pure / denatured		1
Klenk	Ethanol 96% pure / denatured	10059387-15102	1

- 90 688 spectra from 501 *Apo-Ident* customers from a total of 26 645 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Ethanol 96% pure / denatured* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Ethanol 96% pure / denatured* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	160	0	93 614
Type C	0	637	27	90 649

The substance/substance group *Ethanol 96% pure / denatured* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 96.2500 %)
Type C	100.0000 % (> 99.6588 %)	95.9337 % (> 95.4819 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31488	31488	0.00	53.06
31564	31488	4.56	49.97

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Ethyl acetate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30719-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Ethyl acetate; Aether aceticus; Ethylis acetate

Special notes

When selecting the *Ethyl acetate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Ethyl acetate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Ethyl acetate	11357003	30719	40	1402681
Caelo	Ethyl acetate	13366102	31385	40	20131107
Caelo	Ethyl acetate	180804	34075	40	20180320*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Ethyl acetate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 80 spectra from 2 reference samples of the substance/substance group *Ethyl acetate*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Ethyl acetate	15362301	32568	40
Caelo	Ethyl acetate	15362303	33215	40

- 93 794 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 12 spectra from 11 *Apo-Ident* customers from 10 batches from the substance/substance group *Ethyl acetate*.
- These include spectra of independent samples from 9 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Ethyl acetate	15362301	1
Caelo	Ethyl acetate	12357604	1
Caelo	Ethyl acetate	13366101	1
Caelo	Ethyl acetate	15182402	2
Caelo	Ethyl acetate	10286201	1
Caelo	Ethyl acetate	15182404	1
Caelo	Ethyl acetate	9081304	1
Caelo	Ethyl acetate	11357003	2
Caelo	Ethyl acetate	12357601	1
Caelo	Ethyl acetate	15362304	1

- 91 340 spectra from 501 *Apo-Ident* customers from a total of 26 999 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Ethyl acetate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Ethyl acetate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	80	0	93 794
Type C	0	12	0	91 340

The substance/substance group *Ethyl acetate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 92.5000 %)
Type C	100.0000 % (> 99.6593 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to

exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30719	30719	0.00	55.22
31385	34075	7.17	47.97
34075	34075	0.00	56.12

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Eucalyptus oil / Rosmary oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30875-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Eucalyptus oil / Rosmary oil; Eucalypti aetheroleum; Eucalyptus oil; Oleum eucalypti; Oleum rosmarini; Rosemary oil; Rosmarini aetheroleum

Special notes

When selecting the *Eucalyptus oil / Rosmary oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Eucalyptus oil / Rosmary oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Rosemary oil	12E08-N05	30875	40	1502124
Bombastus	Eucalyptus oil	279796	31110	40	73399

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Eucalyptus oil	13106314	31284	34	AR-15-FG-000909-01
Caelo	Eucalyptus oil	14039604	31285	38	1412280
Caelo	Eucalyptus oil	14039609	31298	40	20140214
Euro OTC	Rosemary oil	1303016-01	31311	40	1502125
Caelo	Eucalyptus oil	14121904	31314	40	AR-15-FG-004302-01
Euro OTC	Eucalyptus oil	1406028-01	31527	40	20140716
Euro OTC	Rosemary oil	1511016	31822	40	20151214*
Euro OTC	Eucalyptus oil	1511015	31823	40	20151214*
Caelo	Eucalyptus oil	14039626	31864	40	20140211
Caelo	Eucalyptus oil	15316504	31865	40	20150925
Fagron	Eucalyptus oil	16C15-B05-331849	33277	40	20160411

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 512 spectra of 13 reference samples from the substance/substance group *Eucalyptus oil* / *Rosemary oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 13 different batches.
- 35 032 spectra from a total of 827 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 500 spectra from 15 reference samples of the substance/substance group *Eucalyptus oil* / *Rosemary oil*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Eucalyptus oil	92752399	30253	40
Caelo	Eucalyptus oil	11104312	30859	40
Euro OTC	Eucalyptus oil	1406028-01	31527 [†]	20
Euro OTC	Rosemary oil	1511016	31822 [†]	20
Euro OTC	Eucalyptus oil	1511015	31823 [†]	20
Caelo	Eucalyptus oil	14039626	31864 [†]	20
Caelo	Eucalyptus oil	15316504	31865 [†]	20
Caelo	Rosemary oil	162769	32686	40
Euro OTC	Eucalyptus oil	1608042	32773	40
Euro OTC	Rosemary oil	1608043	32774	40
Caelo	Eucalyptus oil	163513	33084	40
Euro OTC	Eucalyptus oil	1703030	33237	40
Euro OTC	Eucalyptus oil	1801002	33928	40
Euro OTC	Rosemary oil	1803038	34088	40
Euro OTC	Eucalyptus oil	1806016	34189	40

- 93 374 spectra from a total of 2152 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 570 spectra from 185 *Apo-Ident* customers from 262 batches from the substance/substance group *Eucalyptus oil* / *Rosemary oil*.
- These include spectra of independent samples from 248 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Eucalyptus oil	13106326	3
Caelo	Eucalyptus oil	1214E-01978	1
Fagron	Eucalyptus oil	13E15-N04	3
Caelo	Eucalyptus oil	14039602	1
Caelo	Eucalyptus oil	14039609	17
VDL 26.04.14 Preis 5,25	Eucalyptus oil		1
Caelo	Eucalyptus oil	290094	1
Bombastus	Eucalyptus oil	286454	1
Bombastus	Eucalyptus oil	293175	1
Caelo	Eucalyptus oil	14039625	5
Fagron	Eucalyptus oil	14A17-B40	2
Euro OTC	Eucalyptus oil	1211028-02	2
Fagron	Eucalyptus oil	11D11-N02	1
Fagron	Eucalyptus oil	192680-0002	1
Bombastus	Eucalyptus oil	269495	1
Bombastus	Eucalyptus oil	275646	2
Caelo	Eucalyptus oil	276113	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Eucalyptus oil	916M-01978	1
Fagron	Eucalyptus oil	11K23-N05	2
Caelo	Eucalyptus oil	12114112	5
Fagron	Eucalyptus oil	12G25-N11	4
Frey & Lau GmbH/Phoenix	Eucalyptus oil	12G25-N11	1
Caelo	Eucalyptus oil	13106411	3
Caelo	Eucalyptus oil	13106403	3
Caelo	Eucalyptus oil	13106412	3
Caelo	Eucalyptus oil	13106413	1
Fagro/Gehe	Eucalyptus oil	14/20-B02-302706	1
Caelo	Eucalyptus oil	14039601	2
Caelo	Eucalyptus oil	14039603	1
ahd/lamotte	Eucalyptus oil	14039604	1
Caelo	Eucalyptus oil	14039604	1
sana 05.12.14 5,25	Eucalyptus oil	14039609	1
Caelo	Eucalyptus oil	14039610	1
Caelo	Eucalyptus oil	14039618	3
Caelo	Eucalyptus oil	14039622	3
Caelo	Eucalyptus oil	14039626	7
14039625	Eucalyptus oil	14039625	1
Caelo	Eucalyptus oil	14039627	2
Caelo	Eucalyptus oil	14039630	1
Gehe 20.02.15	Eucalyptus oil	1406028-02	1
Euro OTC	Eucalyptus oil	1406028-01	1
Caelo	Eucalyptus oil	14121904	6
Caelo	Eucalyptus oil	1439603	1
Caelo	Eucalyptus oil	14121905	1
Phönix	Eucalyptus oil	14A17-B40	2
Fagron	Eucalyptus oil	14A17-B40-291551	1
Fagron	Eucalyptus oil	14e09-b06-294679	1
Phönix	Eucalyptus oil	14A17-B40-291551	1
phönix	Eucalyptus oil	14E09-B06-294681	1
Caelo	Eucalyptus oil	15011701	5
Fagron	Eucalyptus oil	14E09-B06-294681	1
Caelo	Eucalyptus oil	15011702	1
Caelo	Eucalyptus oil	15011703	1
Caelo	Eucalyptus oil	15011704	1
Caelo	Eucalyptus oil	15063301	4
Caelo	Eucalyptus oil	15063303	1
Caelo	Eucalyptus oil	15063302	1
Caelo	Eucalyptus oil	15063309	13
Caelo	Eucalyptus oil	15063312	9
Euro OTC	Eucalyptus oil	151101501	1
Euro OTC	Eucalyptus oil	1511015-01	2
Caelo	Eucalyptus oil	1515Q-01978	1
Caelo	Eucalyptus oil	15316501	7
Caelo	Eucalyptus oil	15316502	6
Caelo	Eucalyptus oil	15316504	10
Caelo	Eucalyptus oil	15316505	1
Caelo	Eucalyptus oil	15316512	3
Caelo	Eucalyptus oil	15326512	1
Caelo	Eucalyptus oil	15316514	1
Caelo	Eucalyptus oil	15400202	5
Sanacorp	Eucalyptus oil	15400202	1
Celo	Eucalyptus oil	15400202	1
Caelo	Eucalyptus oil	15400203	2
Allianz Healthcare	Eucalyptus oil	15400202	1
Caelo	Eucalyptus oil	15400210	5
Phoenix WE:27.01.2017 ...	Eucalyptus oil	15400210	1
Caelo	Eucalyptus oil	15400211	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Eucalyptus oil	15B03-B02-304362	3
Euro OTC	Eucalyptus oil	1608042-01	2
Caelo	Eucalyptus oil	1608042	2
Ichthyol-Gesellschaft	Eucalyptus oil	1610280E/15316504	2
Caelo	Eucalyptus oil	16235505	3
Caelo	Eucalyptus oil	16351301	2
Sanacorp	Eucalyptus oil	16351301	1
Sanacorp	Eucalyptus oil	16351309	1
Caelo	Eucalyptus oil	16351309	2
Caelo	Eucalyptus oil	16351312	5
ACM	Eucalyptus oil		1
Caelo	Eucalyptus oil	16351313	1
Bombastus	Eucalyptus oil	17000944	1
Bombastus	Eucalyptus oil	17000947	1
Caelo	Eucalyptus oil	17198802	1
Caelo	Eucalyptus oil	28071502	1
Bombastus	Eucalyptus oil	286740	2
Bombastus	Eucalyptus oil	285615	1
Bombastus	Eucalyptus oil	288085	3
Bombastus	Eucalyptus oil	290094	1
Bombastus	Eucalyptus oil	291597	2
Bombastus	Eucalyptus oil	292317	1
Caelo	Eucalyptus oil	292087	1
Bombastus	Eucalyptus oil	293087	2
Bombastus	Eucalyptus oil	293227	3
Caelo	Eucalyptus oil	293087	1
Bombastus	Eucalyptus oil	294181	1
Bombastus	Eucalyptus oil	294460	1
Bombastus	Eucalyptus oil	294687	2
Bombastus	Eucalyptus oil	296079	2
Bombastus	Eucalyptus oil	295573	1
Caelo	Eucalyptus oil	296079	1
Bombastus	Eucalyptus oil	299699	2
Bombastus	Eucalyptus oil	296325	1
Bombastus	Eucalyptus oil	298889	1
Bombastus	Eucalyptus oil	300063	4
Bombastus	Eucalyptus oil	298985	1
Bombastus	Eucalyptus oil	300798	1
Bombastus	Eucalyptus oil	302670	1
Bombastus	Eucalyptus oil	301205	1
Bombastus	Eucalyptus oil	303327	1
L10/Bambastus	Eucalyptus oil	303336	1
Bombastus	Eucalyptus oil	304517	2
Bombastus	Eucalyptus oil	307667	3
Bombastus	Eucalyptus oil	306025	1
Sanacorp	Eucalyptus oil	48/AKD3	1
Taoasis	Eucalyptus oil	8745-115428	1
Caelo	Eucalyptus oil	4136096	2
Sanacorp	Eucalyptus oil		2
Taoasis	Eucalyptus oil	8745-116747	1
Bombastus	Eucalyptus oil	91101	1
Caelo	Eucalyptus oil	513I-01978	1
Caelo	Eucalyptus oil	7011302	1
Caelo	Eucalyptus oil	9351601	1
Caelo	Eucalyptus oil	10051816	1
Caelo	Eucalyptus oil	10318506	1
Fagron	Eucalyptus oil	11102506	1
Euro OTC	Eucalyptus oil	110801001	1
Fagron	Eucalyptus oil	11/27-N03	1
Caesar & Loretz GmbH D-...	Eucalyptus oil	11102507	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
ANZAG 19.12.2012	Eucalyptus oil	110801001	1
Caelo	Eucalyptus oil	11102507	2
Caelo	Eucalyptus oil	11104312	1
Caelo	Eucalyptus oil	11324107	6
Caelo	Eucalyptus oil	11324108	1
Caelo	Eucalyptus oil	11324112	8
Sanacorp	Eucalyptus oil	11324112	1
Euro OTC	Eucalyptus oil	1202030-01	1
Caelo	Eucalyptus oil	12103607	4
Caelo	Eucalyptus oil	12103605	1
Caelo	Eucalyptus oil	12103612	1
Caelo	Eucalyptus oil	12103611	4
Caelo	Eucalyptus oil	12114107	2
Caelo	Eucalyptus oil	12114111	1
Euro OTC	Eucalyptus oil	1211028-01	2
Fagron	Eucalyptus oil	12B23-N08	9
PVL	Eucalyptus oil	1260813	1
Fagron	Eucalyptus oil	12D05-N03	1
Euro OTC	Eucalyptus oil	12D05-N03	1
Caelo	Eucalyptus oil	13106401	4
Spirig/Noweda	Eucalyptus oil	13106401	2
Caelo	Eucalyptus oil	13106301	1
Caelo	Eucalyptus oil	13106402	2
Caelo	Eucalyptus oil	13106304	1
Noveda	Eucalyptus oil	13106402	1
Caelo	Eucalyptus oil	14a17-b40-29155	1
Caelo	Eucalyptus oil	2112I-01978	1
Bombastus	Eucalyptus oil	261706	1
Bombastus	Eucalyptus oil	270004	1
Caelo	Eucalyptus oil	270421	1
Bombastus	Eucalyptus oil	273707	1
Caelo	Eucalyptus oil	275646	1
Anzag 14.06.2013	Eucalyptus oil	275484	1
Anzag	Eucalyptus oil	273707	1
Bombastus	Eucalyptus oil	278619	2
Euro OTC	Eucalyptus oil	278619	2
Bombastus	Eucalyptus oil	276113	1
Bombastus	Eucalyptus oil	277671	1
Bombastus	Eucalyptus oil	279798	1
Bombastus	Eucalyptus oil	281667	2
Bombastu-Werke AG/Gehe	Eucalyptus oil	282165	1
AHD 03.12.2013	Eucalyptus oil	282165	1
Bombastus	Eucalyptus oil	283384	1
Bombastu-Werke AG/Gehe	Eucalyptus oil	283991	1
Bombastus	Eucalyptus oil	285362	2
Caelo	Eucalyptus oil	3012E-01978	1
Bombastus	Eucalyptus oil	284861	1
Taoasis	Eucalyptus oil	502745-108325	1
Taoasis GmbH	Eucalyptus oil	33364-110969	1
Taoasis	Eucalyptus oil	521399-108557	1
Taoasis	Eucalyptus oil	53970-110969	1
Taoasis	Eucalyptus oil	543495-112704	1
Taoasis	Eucalyptus oil	527885-109905	1
Taoasis	Eucalyptus oil	557095-114030	1
Gehe 12.04.2013	Eucalyptus oil		1
Caelo	Eucalyptus oil		1
Bergland / Gehe	Eucalyptus oil	K301874	1
Bergland / Gehe	Eucalyptus oil	K302509	1
Bombastus	Eucalyptus oil		2
Caelo	Eucalyptus oil	13106316	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Eucalyptus oil	14039612	1
VDL;5,40EUR;11.04.15	Eucalyptus oil	14039625	1
sanacorp, 5,03EUR , 5.09....	Eucalyptus oil	15400202	3
Caelo	Eucalyptus oil	15316507	1
Sanacorp, 4,99EUR, 20.1.17	Eucalyptus oil	15400210	2
Bomastus direkt	Eucalyptus oil	293624	1
Fagron	Rosemary oil	12E08-N05	4
Fagron	Rosemary oil	13I19-B04	2
Caelo	Rosemary oil	2613M-02012	1
Euro OTC	Rosemary oil	1202025-01	2
Bombastus	Rosemary oil	296037	1
Bombastus	Rosemary oil	282541	1
Sanacorp, 21.02.17, 8,73EUR	Rosemary oil	15442017	1
Caelo	Rosemary oil	11287208	4
Caelo	Rosemary oil	1209034-01	2
Fagron	Rosemary oil	13B18-N01	3
Caelo	Rosemary oil	12246501	5
Bombastus	Rosemary oil	281945	1
Caelo	Rosemary oil	12246517	2
Euro OTC	Rosemary oil	1209034-01	2
Bergland/Sanacorp	Rosemary oil	C301101-3448311132	1
Bombastus	Rosemary oil	295093	1
Bombastus	Rosemary oil	288951	1
Euro OTC	Rosemary oil	1404037-01	1
Caelo	Rosemary oil	11287211	7
Bombastus	Rosemary oil	273724	3
Caelo	Rosemary oil	11031402	1
Caelo	Rosemary oil	12246520	3
Frey & Lau GmbH/Phoenix	Rosemary oil	13/19-B04-293575	1
Fagron	Rosemary oil	13/19-B04-293577	5
Caelo	Rosemary oil	13/19-B04-293577	3
Euro OTC	Rosemary oil	130801401	2
Fagron	Rosemary oil	13/19-B04-308855	1
Phönix	Rosemary oil	131119-B04	1
Fagron	Rosemary oil	13I19B04	1
Fagron	Rosemary oil	13I19-B04-293577	2
Caelo	Rosemary oil	14092904	1
Euro OTC	Rosemary oil	1404037-02	1
Caelo	Rosemary oil	13I19-B04-293599	1
Caelo	Rosemary oil	14092910	4
Caelo	Rosemary oil	14092911	4
Caelo	Rosemary oil	15442017	5
Caelo	Rosemary oil	15316703	1
Caelo	Rosemary oil	15316704	3
Caelo	Rosemary oil	18091402	1
Taoasis/Noweda	Rosemary oil	1921019-121067	1
Taoasis	Rosemary oil	26057-118014	1
Bombastus	Rosemary oil	285668	1
Bombastus	Rosemary oil	2872889	1
Bombastus	Rosemary oil	278770	1
Bombastus	Rosemary oil	289392	2
Caelo	Rosemary oil	31031502	1
Caelo	Rosemary oil	296036	1
Taoasis/Noweda	Rosemary oil	3231015-118416	1
Taoasis	Rosemary oil	3231031-117639	1
Phoenix/ Taoasis	Rosemary oil	3231016	1
Noweda/Phoenix/ Taoasis	Rosemary oil	3231031-117639	1
Caelo	Rosemary oil	4101301	1
Caelo	Rosemary oil	3091206	1
noweda	Rosemary oil		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PVL	Rosemary oil	491212	1
Noweda	Rosemary oil	10374817	1
Caelo	Rosemary oil	10374817	2
Hedinger	Rosemary oil	11287208	1
Anzag 05.11.12 ek:7,37	Rosemary oil	11287211	1
Noweda	Rosemary oil	11287211	1
Fagron	Rosemary oil	11K21-N05	3
Noweda	Rosemary oil	11287211ui	1
Taoasis	Rosemary oil	121217-110772	1
Caelo	Rosemary oil	12236401	1
Anzag 12.03.13 ek: 7,37EUR	Rosemary oil	12246508	1
Caelo	Rosemary oil	12246508	6
Fagron	Rosemary oil	12246508	1
AHD , 03.09.13	Rosemary oil	12246515	2
Caelo	Rosemary oil	12246515	3
Taoasis GmbH	Rosemary oil	130815-113178	1
Euro OTC	Rosemary oil	1303016-01	1
Gehe 14.06.2013	Rosemary oil	12G12-N12	1
Taoasis	Rosemary oil	130313-111912	1
Fagron	Rosemary oil	13/19-B04-287289	1
Fagron	Rosemary oil	13A28-N03	3
Großhandel	Rosemary oil	13A28-N03	1
fargo	Rosemary oil	13A28-N03	2
Caelo	Rosemary oil	1511016-01	1
Euro OTC	Rosemary oil	1511016-02	1
Caelo	Rosemary oil	15316706	1
Caelo	Rosemary oil	15442001	1
Caelo	Rosemary oil	15442018	1
Caelo	Rosemary oil	16276903	1
Caelo	Rosemary oil	16276904	1
Caelo	Rosemary oil	16276908	2
rosma	Rosemary oil	17000647	1
Caelo	Rosemary oil	17012014	1
Caelo	Rosemary oil	17326702	3
Fagron	Rosemary oil	16L22B84329543	1
Caelo	Rosemary oil	273098	1
Caelo	Rosemary oil	2715E-02012	1
Caelo	Rosemary oil	275963	1
Bombastus	Rosemary oil	278907	4
Caelo	Rosemary oil	279984	1
phönix 14.08.2013	Rosemary oil	275963	1
Caelo	Rosemary oil	28051309	1
Phönix	Rosemary oil	282665	1
Bombastus	Rosemary oil	304472	1
Euro OTC	Rosemary oil	290203901	1
Bombastus	Rosemary oil	301180	1
Bombastus	Rosemary oil	306506	2
Taoasis GmbH	Rosemary oil	3231016-115522	1
Fagron	Rosemary oil	693655-N12	1
Taoasis GmbH	Rosemary oil	7028-105285	1

- 90 782 spectra from 501 *Apo-Ident* customers from a total of 26 751 batches from a further 332 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Eucalyptus oil* / *Rosmary oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Eucalyptus oil* / *Rosmary oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	512	0	35 032
Type B	0	496	4	93 374
Type C	0	556	14	90 782

The substance/substance group *Eucalyptus oil* / *Rosmary oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.8281 %)
Type B	100.0000 % (> 99.9761 %)	99.2000 % (> 98.6000 %)
Type C	100.0000 % (> 99.6586 %)	97.5439 % (> 97.0175 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30875	30875	0.00	41.99
31110	31823	4.82	41.42
31284	31284	0.00	39.39
31285	31285	0.00	40.17
31298	31285	1.48	39.90
31311	31311	0.00	40.66
31314	31314	0.00	40.99
31527	31823	3.14	39.74
31822	31822	0.00	41.48
31823	31823	0.00	42.06
31864	31823	2.35	41.93
31865	31314	5.90	40.86
33277	31823	11.32	39.08

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Eugenol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31286-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Eugenol; Eugenolum

Special notes

When selecting the *Eugenol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Eugenol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Eugenol	13401501	31286	40	1407436
Caelo	Eugenol	13401506	31397	40	20131213
Caelo	Eugenol	173534	33896	40	20180126*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Eugenol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 80 spectra from 2 reference samples of the substance/substance group *Eugenol*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Eugenol	15H11-B03-322946	32731	40
Fagron	Eugenol	15H11-B03-322946	33220	40

- 93 794 spectra from a total of 2166 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 8 spectra from 7 *Apo-Ident* customers from 6 batches from the substance/substance group *Eugenol*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Eugenol	13H14-B40-286617	1
Caelo	Eugenol	14165401	1
Caelo	Eugenol	14165403	1
Caelo	Eugenol	15255901	2
Caelo	Eugenol	16034701	1
Caelo	Eugenol	15255904	2

- 91 344 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Eugenol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Eugenol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	80	0	93 794
Type C	0	8	0	91 344

The substance/substance group *Eugenol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 92.5000 %)
Type C	100.0000 % (> 99.6597 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31286	31286	0.00	37.70
31397	31286	1.56	39.66
33896	33896	0.00	39.23

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Fennel honey (conserved)**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30900-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Fennel honey (conserved); Mel foeniculi (konserviert)

Special notes

When selecting the *Fennel honey (conserved)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Fennel honey (conserved)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Fennel honey (co...	12247203	30900	40	1402574
Caelo	Fennel honey (co...	14045202	31598	40	20140317
Caelo	Fennel honey (co...	181580	34233	40	20180704*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Fennel honey (conserved)*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 80 spectra from 3 reference samples of the substance/substance group *Fennel honey (conserved)*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Fennel honey (conserved)	12247203	30900 [†]	20
Caelo	Fennel honey (conserved)	14045202	31598 [†]	20
Caelo	Fennel honey (conserved)	16039406	33153	40

- 93 794 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2 spectra from 2 *Apo-Ident* customers from 2 batches from the substance/substance group *Fennel honey (conserved)*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Fennel honey (conserved)	13152606	1
Caelo	Fennel honey (conserved)	13152605	1

- 91 350 spectra from 501 *Apo-Ident* customers from a total of 27 007 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Fennel honey (conserved)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Fennel honey (conserved)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	79	1	93 794
Type C	0	2	0	91 350

The substance/substance group *Fennel honey (conserved)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9762 %)	98.7500 % (> 95.0000 %)
Type C	100.0000 % (> 99.6630 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30900	30900	0.00	17.04
31598	30900	10.33	21.23
34233	34233	0.00	21.42

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Fennel oil, bitter**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30398-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Fennel oil, bitter; Oleum foeniculi amari

Special notes

When selecting the *Fennel oil, bitter* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Fennel oil, bitter*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Fennel oil, bitter	14190402	31706	40	20140625
Caelo	Fennel oil, bitter	15160103	31805	40	1703170
Bombastus	Fennel oil, bitter	295236	32922	40	1704221

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Fennel oil, bitter*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 200 spectra from 6 reference samples of the substance/substance group *Fennel oil, bitter*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Fennel oil, bitter	09366708	30169	40
Bombastus	Fennel oil, bitter	264987	30398	40
Caelo	Fennel oil, bitter	14190402	31706 [†]	20
Caelo	Fennel oil, bitter	15160103	31805 [†]	20
Caelo	Fennel oil, bitter	16008606	32961	40
Caelo	Fennel oil, bitter	171809	33680	40

- 93 674 spectra from a total of 2161 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 39 spectra from 25 *Apo-Ident* customers from 21 batches from the substance/substance group *Fennel oil, bitter*.
- These include spectra of independent samples from 20 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Fennel oil, bitter	14190402	2
Caelo	Fennel oil, bitter	15160101	2
Caelo	Fennel oil, bitter	16008601	1
Caelo	Fennel oil, bitter	16008606	1
Caelo	Fennel oil, bitter	12218103	4
Caelo	Fennel oil, bitter	17337601	1
Erika Apotheke	Fennel oil, bitter	12218104	1
Caelo	Fennel oil, bitter	24101403	1
Caelo	Fennel oil, bitter	12218104	5
Caelo	Fennel oil, bitter	12370303	1
Caelo	Fennel oil, bitter	12061412	1
Caelo	Fennel oil, bitter	15160106	1
Caelo	Fennel oil, bitter	12370306	4
Caelo	Fennel oil, bitter	12370312	3
Celo / Kehr	Fennel oil, bitter	12370314	1
Bombastus	Fennel oil, bitter	291598	1
Bombastus	Fennel oil, bitter	278909	3
Caelo	Fennel oil, bitter	268017	1
Sanacorp	Fennel oil, bitter	278909	1
Bombastus	Fennel oil, bitter	285544	1
Bombastus	Fennel oil, bitter	274253	1
Bombastus	Fennel oil, bitter	289232	1
Bombastus	Fennel oil, bitter	281715	1

- 91 313 spectra from 501 *Apo-Ident* customers from a total of 26 988 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Fennel oil, bitter* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Fennel oil, bitter* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	170	30	93 674
Type C	0	35	4	91 313

The substance/substance group *Fennel oil, bitter* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	85.0000 % (> 83.5000 %)
Type C	100.0000 % (> 99.6588 %)	89.7436 % (> 82.0513 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31706	31805	1.54	62.45
31805	31805	0.00	62.22
32922	32922	0.00	45.64

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Gel Cordes®
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30632-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Gel Cordes®

Special notes

When selecting the *Gel Cordes®* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Gel Cordes®*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Ichthyol	Gel Cordes®	12A035	30632	40	1405096
Ichthyol	Gel Cordes®	13A007	31083	40	20130513
Ichthyol	Gel Cordes®	15A099	31876	40	20150624
Ichthyol	Gel Cordes®	15A099	32911	40	20150624

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Gel Cordes*[®]. These samples are listed above in the section *calibration samples*. The reference samples originate from 3 different batches.
- 35 384 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 176 spectra from 6 reference samples of the substance/substance group *Gel Cordes*[®].
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Ichthyol	Gel Cordes [®]	14A017	30714	40
Ichthyol	Gel Cordes [®]	15A099	31876 [†]	20
Ichthyol	Gel Cordes [®]	15A099	32910	40
Ichthyol	Gel Cordes [®]	17A047	33987	40
Ichthyol	Gel Cordes [®]	15A099	not required	24
Ichthyol	Gel Cordes [®]	16A169	not required	12

- 93 698 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 24 spectra from 12 *Apo-Ident* customers from 6 batches from the substance/substance group *Gel Cordes*[®].
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Ichthyol/Anzag	Gel Cordes [®]	12A035	2
Caelo	Gel Cordes [®]	13A007	2
Ichthyol-Gesellschaft C...	Gel Cordes [®]	13A007	1
Noweda/Ichthyol-Gesells...	Gel Cordes [®]	14A017	3
Ichthyol/Sanacorp	Gel Cordes [®]	14A017	1
Noweda/Ichthyol-Gesells...	Gel Cordes [®]	15A099	1
Ichthyol/Gehe	Gel Cordes [®]	15A099	1
Ichthyol/Alliance Healt...	Gel Cordes [®]	15A099	3
Ichthyol-Gesellschaft/N...	Gel Cordes [®]	15A099	2
Ichthyol	Gel Cordes [®]	15A099	1
Ichthyol/Sanacorp	Gel Cordes [®]	16A169	1
Ichthyol Gesellschaft	Gel Cordes [®]	17A047	1
AHD	Gel Cordes [®]	16A169	1
Ichthyol	Gel Cordes [®]	17A047	2
Ichthyol/Sanacorp	Gel Cordes [®]	17A047	1
Ichthyol/Phönix	Gel Cordes [®]	17A047	1

- 91 328 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Gel Cordes*[®] can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Gel Cordes*[®] and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	136	40	93 698
Type C	0	24	0	91 328

The substance/substance group *Gel Cordes*[®] can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	77.2727 % (> 75.5682 %)
Type C	100.0000 % (> 99.6589 %)	100.0000 % (> 75.0000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30632	30632	0.00	40.76
31083	30632	15.11	40.10
31876	30632	15.51	40.29
32911	30632	12.10	40.35

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Geranium oil, genuine
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30468-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Geranium oil, genuine; Oleum geranii verum

Special notes

When selecting the *Geranium oil, genuine* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Geranium oil, genuine*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Geranium oil, ge...	12354109	31193	40	1403533
Caelo	Geranium oil, ge...	15172205	31789	40	20150608

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Geranium oil, genuine*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Geranium oil, genuine*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Geranium oil, genuine	11275601	30468	40
Caelo	Geranium oil, genuine	15172205	31789 [†]	20
Caelo	Geranium oil, genuine	171503	33316	40
Caelo	Geranium oil, genuine	16070404	33438	40
Caelo	Geranium oil, genuine	180259	34053	40

- 93 694 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 4 spectra from 4 *Apo-Ident* customers from 4 batches from the substance/substance group *Geranium oil, genuine*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Geranium oil, genuine	16070403	1
Caelo	Geranium oil, genuine	11275605	1
Caelo	Geranium oil, genuine	153670202	1
Caelo	Geranium oil, genuine	35434343	1

- 91 348 spectra from 501 *Apo-Ident* customers from a total of 27 005 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Geranium oil, genuine* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Geranium oil, genuine* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	180	0	93 694
Type C	0	4	0	91 348

The substance/substance group *Geranium oil, genuine* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.6667 %)
Type C	100.0000 % (> 99.6608 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31193	31193	0.00	75.66
31789	31193	4.72	73.58

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Glucose syrup**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30934-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Glucose syrup; Glucose syrup GE 40%; Glucosum liquidum; Sirupus Glucosi

Special notes

When selecting the *Glucose syrup* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Glucose syrup*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Glucose syrup	12I18-N11	30934	40	1403229
Fagron	Glucose syrup	14G16-T01-018036	31603	40	20150120
Fagron	Glucose syrup	16I02-T07-039136	33632	40	20161222

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Glucose syrup*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Glucose syrup*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Glucose syrup	14G16-T01-018036	31603 [†]	20
Fagron	Glucose syrup	14G16-T01-018036	32272	40
Fagron	Glucose syrup	16B02-T11-029595	32273	40

- 93 774 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 9 spectra from 7 *Apo-Ident* customers from 6 batches from the substance/substance group *Glucose syrup*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Glucose syrup	14G16-T01-018036	1
Fagron	Glucose syrup	10J25-T52-013485	1
Fagron	Glucose syrup	16I02-T07-039136	2
Fagron	Glucose syrup	16B02-T11-029595	1
Fagron	Glucose syrup	13D23-N02	1
Fagron	Glucose syrup	13G03-N01	2
Phönix	Glucose syrup	13D23-N02	1

- 91 343 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Glucose syrup* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Glucose syrup* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	100	0	93 774
Type C	0	9	0	91 343

The substance/substance group *Glucose syrup* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6595 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30934	30934	0.00	37.07
31603	30934	6.57	35.88
33632	30934	3.00	33.75

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Glycerine ointment (conserved)**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30986-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Glycerine ointment (conserved); Glyceroli unguentum (konserviert); Unguentum glycerini (konserviert)

Special notes

When selecting the *Glycerine ointment (conserved)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Glycerine ointment (conserved)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Glycerine ointme...	12288201	30986	40	20121017
Caelo	Glycerine ointme...	14031701	31296	40	1410452

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Glycerine ointment (conserved)*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Glycerine ointment (conserved)*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Glycerine ointment (conserved)	16192001	32589	40
Caelo	Glycerine ointment (conserved)	17049101	33240	40
Caelo	Glycerine ointment (conserved)	18095401	34206	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 15 spectra from 6 *Apo-Ident* customers from 7 batches from the substance/substance group *Glycerine ointment (conserved)*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Glycerine ointment (conserved)	13250102	2
Caelo	Glycerine ointment (conserved)	14183401	6
Caelo	Glycerine ointment (conserved)	16192001	2
Caelo	Glycerine ointment (conserved)	15025101	2
Caelo	Glycerine ointment (conserved)	16021501	1
Caelo	Glycerine ointment (conserved)	16192002	1
Caelo	Glycerine ointment (conserved)	14031701	1

- 91 337 spectra from 501 *Apo-Ident* customers from a total of 27 002 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Glycerine ointment (conserved)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Glycerine ointment (conserved)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	120	0	93 754
Type C	0	15	0	91 337

The substance/substance group *Glycerine ointment (conserved)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6592 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30986	31296	11.02	34.97
31296	31296	0.00	35.65

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Glycerol >99%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30189-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Glycerol >99%; Glycerinum anhydricum >99%

Special notes

When selecting the *Glycerol >99%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Glycerol >99%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Glycerol >99%	12111504	30871	40	20120418
Caelo	Glycerol >99%	13127603	31152	40	1405278
Caelo	Glycerol >99%	13370502	31352	40	20131115
Caelo	Glycerol >99%	162308	32693	40	20160809*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Glycerol* >99%. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 309 spectra from 7 reference samples of the substance/substance group *Glycerol* >99%.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Glycerol >99%	10198706	30189	89
Caelo	Glycerol >99%	10198706	30231	40
Caelo	Glycerol >99%	13370502	31352 [†]	20
Caelo	Glycerol >99%	162308	32325	40
Caelo	Glycerol >99%	171857	33703	40
Caelo	Glycerol >99%	171857	33708	40
Caelo	Glycerol >99%	180984	34064	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 565 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 413 spectra from 153 *Apo-Ident* customers from 66 batches from the substance/substance group *Glycerol >99%*.
- These include spectra of independent samples from 63 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Glycerol >99%	13127610	3
Caelo	Glycerol >99%	13370502	9
Caelo	Glycerol >99%	13370503	15
Caelo	Glycerol >99%	14373403	19
Caelo	Glycerol >99%	14373409	15
Caelo	Glycerol >99%	14373413	27
Caelo	Glycerol >99%	1474303	1
Beiersdorf AG	Glycerol >99%	15406501	1
Caelo	Glycerol >99%	15406501	13
Caelo	Glycerol >99%	16023506	21
Caelo	Glycerol >99%	16023504	18
Caelo	Glycerol >99%	17185703	4
Caelo	Glycerol >99%	11217403	1
Caelo	Glycerol >99%	11305101	1
Caelo	Glycerol >99%	13002902	3
Caesar & Lorenz	Glycerol >99%	13002904	1
Caelo	Glycerol >99%	733514373412	1
Caelo	Glycerol >99%	4031501	1
Caelo	Glycerol >99%	13127604	10
Caelo	Glycerol >99%	13127607	2
Caelo	Glycerol >99%	13111403	1
Caelo	Glycerol >99%	13127608	6
Caelo	Glycerol >99%	13370501	2
Caelo	Glycerol >99%	133370507	1
Phönix	Glycerol >99%	13127610	1
Gehe 16.12.14	Glycerol >99%	13370503	1
Caelo	Glycerol >99%	14074302	17
Caelo	Glycerol >99%	13370507	1
Caesar & Loretz GmbH Hi...	Glycerol >99%	14074302	1
Fagron	Glycerol >99%	14074302	1
Caelo	Glycerol >99%	14074303	4
Caelo	Glycerol >99%	14373401	2
Caesar & Loretz GmbH	Glycerol >99%	14373403	1
Caelo	Glycerol >99%	14373404	3
Noweda	Glycerol >99%	14373403	1
Fagron	Glycerol >99%	14373409	2
Caelo	Glycerol >99%	14373412	3
Noweda	Glycerol >99%	14373413	1
Caelo	Glycerol >99%	1473403	1
Caelo	Glycerol >99%	1473409	1
Caelo	Glycerol >99%	15283603	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol >99%	15406502	17
Caesar & Loretz GmbH	Glycerol >99%	15406502	1
Caelo	Glycerol >99%	15406503	15
Caelo	Glycerol >99%	15406504	10
Caelo	Glycerol >99%	16023503	7
Euro OTC	Glycerol >99%	16023503	1
Phoenix	Glycerol >99%	16023504	2
Caesar & Loretz GmbH Hi...	Glycerol >99%	16023506	1
GEHE	Glycerol >99%	16023506	1
Caelo	Glycerol >99%	16023508	5
Caelo	Glycerol >99%	16230802	4
Caelo	Glycerol >99%	16230804	14
Caesar & Loretz GmbH	Glycerol >99%	16230804	1
Caelo	Glycerol >99%	16230805	1
Caelo	Glycerol >99%	16320804	1
Caelo	Glycerol >99%	17037503	2
Caelo	Glycerol >99%	17037504	16
Caelo	Glycerol >99%	17041403	1
Caelo	Glycerol >99%	17185702	1
Caelo	Glycerol >99%	17185704	1
Caelo	Glycerol >99%	18081504	1
Caelo	Glycerol >99%	17292803	1
Caelo	Glycerol >99%	2061201	1
Caelo	Glycerol >99%	2114Q1590	1
Caelo	Glycerol >99%	11061303	1
Caelo	Glycerol >99%	1004211590	1
Caelo	Glycerol >99%	510211590	1
Caelo	Glycerol >99%	11118104	1
Caelo	Glycerol >99%	11305102	3
Phönix 22.08.12	Glycerol >99%	11305102	1
Caelo	Glycerol >99%	12014304	15
Gehe	Glycerol >99%	12014304	1
Caelo	Glycerol >99%	12111503	5
Caelo	Glycerol >99%	12111504	3
Gehe	Glycerol >99%	12111504	2
Caelo	Glycerol >99%	12111505	1
Caelo	Glycerol >99%	1211503	1
Caelo	Glycerol >99%	13002903	9
Caesar/Lorenz	Glycerol >99%	13002904	2
Caelo	Glycerol >99%	13002904	5
Phönix	Glycerol >99%	13002904	2
Fagron	Glycerol >99%	13002904	1
Gehe	Glycerol >99%	13002904	1
Caelo	Glycerol >99%	1300291	1
Caelo	Glycerol >99%	13127603	18
Phönix	Glycerol >99%	13127604	1
Caelo	Glycerol >99%	2713M1590	1
Caelo	Glycerol >99%	13002907	1
Gehe	Glycerol >99%	13127603	2
Caesar & Loretz GmbH Hi...	Glycerol >99%	13127608	2
apomix/ Phoenix	Glycerol >99%	13127608	2
Anzag	Glycerol >99%	13127610	2
Caelo	Glycerol >99%	1317604	1
Phönix/EuRho	Glycerol >99%	12111503	2

- 90 939 spectra from 501 *Apo-Ident* customers from a total of 26 945 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are docu-

mented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Glycerol >99%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Glycerol >99%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	309	0	93 565
Type C	0	410	3	90 939

The substance/substance group *Glycerol >99%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.0583 %)
Type C	100.0000 % (> 99.6586 %)	99.2736 % (> 98.5472 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30871	31152	4.08	108.22
31152	31152	0.00	110.88
31352	31152	4.26	103.87
32693	32693	0.00	176.86

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Glycerol 85%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30886-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Glycerol 85%; Glycerinum 85%; Glycerolum 85%

Special notes

When selecting the *Glycerol 85%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Glycerol 85%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Glycerol 85%	12216606	30886	70	20120731
Fagron	Glycerol 85%	12H16-N04	30980	80	20120905
Caelo	Glycerol 85%	13356404	31360	40	20131112
Caelo	Glycerol 85%	13406202	31377	40	20131219
Euro OTC	Glycerol 85%	1409035-01	31588	40	AR-15-FG-007030-02

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Euro OTC	Glycerol 85%	1506036	31719	40	0000010478
Euro OTC	Glycerol 85%	1509023	31737	40	20151005*
Euro OTC	Glycerol 85%	1601059	31956	40	20160225*
Caelo	Glycerol 85%	17197706	34031	40	AR-18-FG-009772-02

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 430 spectra of 9 reference samples from the substance/substance group *Glycerol 85%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 9 different batches.
- 35 114 spectra from a total of 831 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 870 spectra from 25 reference samples of the substance/substance group *Glycerol 85%*.
- These include spectra of independent samples from 16 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Glycerol 85%	12216606	30886 [†]	10
Caelo	Glycerol 85%	13356404	31360 [†]	20
Euro OTC	Glycerol 85%	1409035-01	31588 [†]	20
Euro OTC	Glycerol 85%	1506036	31719 [†]	20
Euro OTC	Glycerol 85%	1509023	31737 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Glycerol 85%	1601059	31956 [†]	20
Caelo	Glycerol 85%	160637	32203	40
Caelo	Glycerol 85%	160637	32314	40
Caelo	Glycerol 85%	162153	32326	40
Caelo	Glycerol 85%	162189	32328	40
Caelo	Glycerol 85%	162978	32694	40
Euro OTC	Glycerol 85%	1607037	32703	40
Euro OTC	Glycerol 85%	1610032	32828	40
Caelo	Glycerol 85%	163550	33097	40
Caelo	Glycerol 85%	163472	33098	40
Caelo	Glycerol 85%	171577	33320	40
Euro OTC	Glycerol 85%	1703003	33449	40
Euro OTC	Glycerol 85%	1708022	33661	40
Caelo	Glycerol 85%	172155	33704	40
Caelo	Glycerol 85%	171977	33705	40
Caelo	Glycerol 85%	171977	33706	40
Caelo	Glycerol 85%	171577	33709	40
Euro OTC	Glycerol 85%	1801025	33996	40
Euro OTC	Glycerol 85%	1804017	34091	40
Euro OTC	Glycerol 85%	1808009	34218	40

- 93 004 spectra from a total of 2145 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2155 spectra from 392 *Apo-Ident* customers from 487 batches from the substance/substance group *Glycerol 85%*.
- These include spectra of independent samples from 479 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	11217302	2
Caelo	Glycerol 85%	11324201	3
Hedinger	Glycerol 85%	40519	3
Caelo	Glycerol 85%	12108201	23
Caelo	Glycerol 85%	465911	1
Noweda/Hedinger	Glycerol 85%	40509	1
Hedinger/Phönix	Glycerol 85%	40513	1
Gehe-Hedinger	Glycerol 85%	40522	1
Hedinger / Gehe	Glycerol 85%	40515	1
Hedinger/Phoenix	Glycerol 85%	40515	1
Hedinger/Noweda	Glycerol 85%	40515	2

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix/Hedinger	Glycerol 85%	40524	1
Lamotte//AHD	Glycerol 85%	125545	1
Euro OTC	Glycerol 85%	40524	1
Hedinger	Glycerol 85%	40531	2
Hedinger	Glycerol 85%	40524	1
Hedinger/Gehe	Glycerol 85%	40529	2
Hedinger	Glycerol 85%	40538	3
Stumpf/Hedinger	Glycerol 85%	40531	1
Noweda/Hedinger	Glycerol 85%	40538	2
Hedinger/Gehe	Glycerol 85%	40540	2
Caelo	Glycerol 85%	40542	3
Hedinger/Phoenix	Glycerol 85%	40542	1
Hedinger/PX	Glycerol 85%	40547	1
Hedinger	Glycerol 85%	40548	2
Hedinger	Glycerol 85%	40556	2
Hedinger/Sanacorp	Glycerol 85%	40548	2
Hedinger/Gehe	Glycerol 85%	40557	1
Hedinger/Gehe	Glycerol 85%	40548	1
Caelo	Glycerol 85%	40556	1
August Hedinger	Glycerol 85%	40557	1
Hedinger	Glycerol 85%	40557	3
AHD / Hedinger	Glycerol 85%	40557	1
Hedinger/Gehe	Glycerol 85%	40559	1
Stumpf / Hedinger	Glycerol 85%	40557	1
Hedinger/Sanacorp	Glycerol 85%	40559	1
Hedinger	Glycerol 85%	40566	1
Hedinger/Phönix	Glycerol 85%	40567	1
Hedinger/ Sanacorp	Glycerol 85%	40567	1
Hedinger/Phönix	Glycerol 85%	40559	1
Noweda/Hedinger	Glycerol 85%	40574	2
Hedinger/Noweda	Glycerol 85%	40567	1
Hedinger/Noweda	Glycerol 85%	40574	1
Hedinger/ Noweda	Glycerol 85%	40574	1
Hedinger/Gehe	Glycerol 85%	40576	1
Hedinger/Phönix	Glycerol 85%	40576	1
Rondell Apotheke	Glycerol 85%	40576	1
Heedinger/Noweda	Glycerol 85%	40576	1
Hedinger/Noweda	Glycerol 85%	40579	2
Hedinger	Glycerol 85%	40576	4
Hedinger/Phoenix	Glycerol 85%	40579	1
Hedinger/Gehe	Glycerol 85%	40584	3
Hedinger/Sanacorp	Glycerol 85%	40579	1
Hedinger	Glycerol 85%	40586	1
Hedinger/ Sanacorp	Glycerol 85%	40584	2
Phoenix	Glycerol 85%	40596	1
Hedinger	Glycerol 85%	40596	1
Hedinger/Phoenix	Glycerol 85%	40596	1
Hedinger	Glycerol 85%	40598	2
Hedinger/ Sanacorp	Glycerol 85%	40599	2
Hedinger	Glycerol 85%	40599	5
Jenne	Glycerol 85%	40599	1
Alte Apotheke	Glycerol 85%	40599	2
Hedinger/Gehe	Glycerol 85%	40599	1
Caelo	Glycerol 85%	40599	1
Gehe/Hedinger	Glycerol 85%	40599	2
Caelo	Glycerol 85%	40607	1
Hedinger	Glycerol 85%	40612	3
Hedinger	Glycerol 85%	40607	3
Hedinger/Noweda	Glycerol 85%	40612	2
Hedinger/Noweda	Glycerol 85%	40607	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Glycerol 85%	40618	2
Hedinger/ Phoenix	Glycerol 85%	40628	1
Hedinger	Glycerol 85%	40620	1
Gehe/Hedinger	Glycerol 85%	40620	1
Phönix	Glycerol 85%	40628	1
Hedinger / Phönix	Glycerol 85%	40628	2
Caelo	Glycerol 85%	40629	2
Caelo	Glycerol 85%	5051501	1
Fischer	Glycerol 85%	1075012	1
Caelo	Glycerol 85%	8121402	1
otto fischar	Glycerol 85%	1075111	1
Hedinger/Sanacorp	Glycerol 85%	61016B	1
Fischar	Glycerol 85%	1075061	1
Fischar/ Gehe	Glycerol 85%	1076011	1
Sanacorp/Fischar	Glycerol 85%	1076021	1
Otto Fischar GmbH & Co...	Glycerol 85%	1077011	1
Fischer/ Noweda	Glycerol 85%	1076111	1
Caelo	Glycerol 85%	12216601	8
Caelo	Glycerol 85%	11355301	1
Caelo	Glycerol 85%	13144604	7
Caelo	Glycerol 85%	13144606	6
Caelo	Glycerol 85%	13205003	1
Caelo	Glycerol 85%	13276024	1
Caelo	Glycerol 85%	13248201	15
noweda	Glycerol 85%	13276202	1
Caelo	Glycerol 85%	13276202	1
Caelo	Glycerol 85%	13276203	18
Caelo	Glycerol 85%	13276204	13
Caelo	Glycerol 85%	13276205	23
noweda	Glycerol 85%	13356403	1
Caelo	Glycerol 85%	13356403	3
Caelo	Glycerol 85%	13356404	13
Caelo	Glycerol 85%	13356405	13
Gehe	Glycerol 85%	13356405	1
Phönix	Glycerol 85%	13356406	1
Caelo	Glycerol 85%	13356406	10
Caelo	Glycerol 85%	13406202	8
Caelo	Glycerol 85%	13406203	1
Caelo	Glycerol 85%	13406204	7
Caelo	Glycerol 85%	13406205	4
Fagron	Glycerol 85%	13I12-B03-285611	1
Fagron	Glycerol 85%	13J01-B01-285781	2
Fagron	Glycerol 85%	13J20-B01-285782	1
Fagron	Glycerol 85%	13J29-B01-285782	1
Euro OTC	Glycerol 85%	1401009-01	8
Euro OTC	Glycerol 85%	1403006-01	8
Euro OTC	Glycerol 85%	14050021-02	1
Euro OTC	Glycerol 85%	1405021-01	4
Euro OTC	Glycerol 85%	1405021-02	7
Euro OTC	Glycerol 85%	1409035-01	18
Caelo	Glycerol 85%	1409035-01	2
Caelo	Glycerol 85%	14148103	5
Caelo	Glycerol 85%	14148105	2
Caelo	Glycerol 85%	14148106	1
Ph.Eur.7.0	Glycerol 85%	14148106	1
Caelo	Glycerol 85%	14149403	19
Caelo	Glycerol 85%	14149404	10
Phönix	Glycerol 85%	14149403	1
Caelo	Glycerol 85%	14149406	12
Caelo	Glycerol 85%	14149410	17

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	14149411	11
Sanacorp	Glycerol 85%	14149410	2
Caelo	Glycerol 85%	1419411	1
Caelo	Glycerol 85%	14194911	1
Caelo	Glycerol 85%	14249301	4
Caelo	Glycerol 85%	1418103	1
Fagron	Glycerol 85%	14306202	2
Caelo	Glycerol 85%	14306202	1
Caelo	Glycerol 85%	14306203	13
Caelo	Glycerol 85%	14306204	9
Caelo	Glycerol 85%	14306205	8
Euro OTC	Glycerol 85%	1502010-01	15
Caelo	Glycerol 85%	1502010-01	1
Caelo	Glycerol 85%	15036902	11
Bombastus	Glycerol 85%	1502010-01	1
Caelo	Glycerol 85%	15036901	1
Caelo	Glycerol 85%	15036903	11
Sanacorp	Glycerol 85%	15036902	1
Phoenix	Glycerol 85%	15036903	1
Fagron	Glycerol 85%	15036903	1
Sanacop	Glycerol 85%	15036903	1
Caelo	Glycerol 85%	15036904	13
Caelo	Glycerol 85%	15050101	2
Caelo	Glycerol 85%	15049801	2
Caelo	Glycerol 85%	15036905	2
Caelo	Glycerol 85%	15050102	6
Euro OTC	Glycerol 85%	1506036-01	16
Euro OTC	Glycerol 85%	1509023-01	11
Caelo	Glycerol 85%	1514604	1
Caelo	Glycerol 85%	15146301	7
Bombastus	Glycerol 85%	15146303	1
Caelo	Glycerol 85%	15115005	2
Caelo	Glycerol 85%	15146304	2
Caelo	Glycerol 85%	15146305	4
Caelo	Glycerol 85%	15283505	1
Caelo	Glycerol 85%	15283503	1
Caelo	Glycerol 85%	15283506	3
Caelo	Glycerol 85%	1528503	1
Sanacorp 09.12.15	Glycerol 85%	152902202	1
Caelo	Glycerol 85%	15283504	1
Caelo	Glycerol 85%	1529201	1
Caelo	Glycerol 85%	15292201	14
Caelo	Glycerol 85%	15292202	11
Fagron	Glycerol 85%	15292201	1
Caelo	Glycerol 85%	15292203	8
Caelo	Glycerol 85%	15292204	20
Sanacop	Glycerol 85%	15292204	1
Henry Lamotte	Glycerol 85%	15292204	1
Caelo	Glycerol 85%	15292205	18
Sanacorp	Glycerol 85%	15292205	2
Caelo	Glycerol 85%	153782021Flasche	1
Caelo	Glycerol 85%	15378202	8
Caelo	Glycerol 85%	153782032Flasche	1
Caelo	Glycerol 85%	15378203	10
Phönix	Glycerol 85%	15378203	2
Caelo	Glycerol 85%	15378210	1
Caelo	Glycerol 85%	15378209	2
Caelo	Glycerol 85%	15379201	4
Caelo	Glycerol 85%	15404502	1
Caelo	Glycerol 85%	15379208	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	15404504	5
Caelo	Glycerol 85%	15404505	4
Euro OTC	Glycerol 85%	1601059-01	27
Caelo	Glycerol 85%	16063701	15
Caelo	Glycerol 85%	16063702	13
Caesar & Loretz GmbH Hi...	Glycerol 85%	16063702	1
Caelo	Glycerol 85%	16063707	2
Caelo	Glycerol 85%	16063706	2
Euro OTC	Glycerol 85%	1607037-01	9
Caelo	Glycerol 85%	1607037-01	1
Euro OTC	Glycerol 85%	1610032-01	12
Caelo	Glycerol 85%	16215301	9
EuRho / Phoenix	Glycerol 85%	1610032-01	1
Caelo	Glycerol 85%	16215302	9
Caelo	Glycerol 85%	16215303	16
Caelo	Glycerol 85%	16215310	3
Caelo	Glycerol 85%	16218901	4
Cordes/Noweda	Glycerol 85%	1627803	1
Caelo	Glycerol 85%	1625303	1
Caelo	Glycerol 85%	16218902	2
Caesar&Loretz/Noweda	Glycerol 85%	16218906	2
Caelo	Glycerol 85%	1629706	1
Caelo	Glycerol 85%	16297802	1
Caelo	Glycerol 85%	16297803	7
Caeolo/Gehe	Glycerol 85%	16297802	3
Caelo	Glycerol 85%	16297806	25
Euro OTC	Glycerol 85%	16297806	1
Caelo	Glycerol 85%	16347202	3
Caelo	Glycerol 85%	16297808	1
Caelo	Glycerol 85%	16347203	2
Caelo	Glycerol 85%	16355001	2
Caelo	Glycerol 85%	16347209	1
Caelo	Glycerol 85%	16355002	4
Caleo/Phönix	Glycerol 85%	16355002	2
Euro OTC	Glycerol 85%	1703003-01	13
Caelo	Glycerol 85%	1703003-01	2
Caelo	Glycerol 85%	17157702	9
Caelo	Glycerol 85%	17061404	1
Bombastus	Glycerol 85%	17157703	1
Caelo	Glycerol 85%	1715772	1
Caelo	Glycerol 85%	17215501	5
Caelo	Glycerol 85%	17157911	1
Caelo	Glycerol 85%	17215502	1
Caelo	Glycerol 85%	17322101	2
Euro OTC	Glycerol 85%	1801025-01	1
Caelo	Glycerol 85%	17322104	1
Caelo	Glycerol 85%	2714M-01589	1
Bombastus	Glycerol 85%	284400	6
Gehe	Glycerol 85%	287120	2
Bombastus	Glycerol 85%	285366	3
Bombastus	Glycerol 85%	287890	2
Bombastus	Glycerol 85%	289295	7
Bombastus	Glycerol 85%	291085	3
Bombastus	Glycerol 85%	292791	17
Bombastus	Glycerol 85%	291086	5
Bombastus	Glycerol 85%	293387	4
Eigenherstellung	Glycerol 85%	295501	1
2959	Glycerol 85%	2959	1
Bombastus	Glycerol 85%	295501	6
Bombastus	Glycerol 85%	294927	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Glycerol 85%	296481	5
Bombastus	Glycerol 85%	297034	1
Bombastus	Glycerol 85%	298592	3
Bombastus	Glycerol 85%	298593	6
Bombastus	Glycerol 85%	300393	1
Bombastus	Glycerol 85%	300394	5
Bombastus	Glycerol 85%	300395	5
Fagron	Glycerol 85%	301816	2
Bombastus	Glycerol 85%	301818	2
Bombastus	Glycerol 85%	305199	1
Bombastus	Glycerol 85%	305200	7
Bombastus	Glycerol 85%	306829	2
Bombastus	Glycerol 85%	310592	2
Fagron	Glycerol 85%	3111	1
Nowed 18.12.2014	Glycerol 85%		2
Fischer	Glycerol 85%	5063091	1
Fischar	Glycerol 85%	5055031	1
Fischar	Glycerol 85%	5063091	1
Fischar	Glycerol 85%	5063092	1
Fischar	Glycerol 85%	5064021	1
Caelo	Glycerol 85%	5015Q-01589	1
Caelo	Glycerol 85%	5063111	1
Fischar/V.d.L.	Glycerol 85%	5064071	1
fischar gmbh	Glycerol 85%	5064011	1
fischar gmbh	Glycerol 85%	5064021	1
Noweda/Ichthyol	Glycerol 85%	5064071	1
Caelo	Glycerol 85%	5064091	1
Fischer/Noweda	Glycerol 85%	5065011	1
Fischar/ Now	Glycerol 85%	5065011	1
Fischar	Glycerol 85%	5065011	1
Fischar / Phönix	Glycerol 85%	5065031	2
Fischar	Glycerol 85%	5065031	2
Caelo	Glycerol 85%	5065111	3
Fischar/ Now	Glycerol 85%	5065111	2
Fischer	Glycerol 85%	5065111	1
Fischar/AHD	Glycerol 85%	5065111	1
Fischer	Glycerol 85%	5066061	2
Fischar	Glycerol 85%	5066021	1
Fischer	Glycerol 85%	50660631	1
AHD/Fischar	Glycerol 85%	5067091	1
Fischar	Glycerol 85%	5066081	2
fischar/gehe	Glycerol 85%	5067111	2
fischar/Noweda	Glycerol 85%	5067011	1
Caelo	Glycerol 85%	5068012	1
Herbeta Arzneimittel Be...	Glycerol 85%		1
Lamotte	Glycerol 85%	5538401	7
Jenne/Lamotte	Glycerol 85%	5603001	1
Lamotte	Glycerol 85%	5603001	7
Henry Lamotte	Glycerol 85%	5603001	1
Lamotte	Glycerol 85%	5603002	1
Fagron	Glycerol 85%	5763701	1
Hanry Lanotte/Phönix	Glycerol 85%	8001116001	1
Sanacorp	Glycerol 85%	79/AKD12	1
Sanacorp	Glycerol 85%	78/AKD7	1
Henry Lamotte	Glycerol 85%	8002208001	1
Noweda	Glycerol 85%	272305	1
spinnrad	Glycerol 85%	914409013	1
Test	Glycerol 85%	Test	1
Caelo	Glycerol 85%	271511	1
Bombastus	Glycerol 85%	272305	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Glycerol 85%	273719	1
Caelo	Glycerol 85%	123142	2
Henry Lamotte	Glycerol 85%	123142	2
Lamotte	Glycerol 85%	123142	1
Lamotte	Glycerol 85%	123304	4
Henry Lamotte	Glycerol 85%	123304	1
Henry	Glycerol 85%	123304	2
Lamotte	Glycerol 85%	123464	5
Sanacorp/ Hedding	Glycerol 85%	400473	2
Hedinger Jenne	Glycerol 85%	40462	1
Hedinger	Glycerol 85%	40464	2
Hedinger Phönix	Glycerol 85%	40471	1
Hedinger	Glycerol 85%	40473	1
Hedinger/Gehe	Glycerol 85%	40473	1
Hedinger	Glycerol 85%	40471	2
Großhandel	Glycerol 85%	40473	1
Anzag/Seeger	Glycerol 85%	40475	2
Sanacorp	Glycerol 85%	40480	1
Hedinger/ Sanacorp	Glycerol 85%	40479	2
Anzag/ Hedinger	Glycerol 85%	40478	1
Caelo	Glycerol 85%	40481	1
Hedinger/Anzag	Glycerol 85%	40483	2
Hedinger	Glycerol 85%	40483	2
Caelo	Glycerol 85%	40483	1
Hedinger/ Anzag	Glycerol 85%	40491	1
Phoenix	Glycerol 85%	40495	2
Sanacorp	Glycerol 85%	40495	1
Hedinger	Glycerol 85%	40497	5
Phönix/Hedinger	Glycerol 85%	40497	1
Sanacorp	Glycerol 85%	40506	1
Hedinger/Phönix	Glycerol 85%	40501	1
Caelo	Glycerol 85%	40507	1
Hedinger	Glycerol 85%	40506	2
Hedinger	Glycerol 85%	40507	1
Hedinger	Glycerol 85%	40509	1
Bartel/Sanacorp	Glycerol 85%	40515	1
Caelo	Glycerol 85%	1113E-01589	1
Caelo	Glycerol 85%	613Q-01589	1
Caelo	Glycerol 85%	607211589	1
Europharm/Otto	Glycerol 85%	1108032-01	1
Caelo	Glycerol 85%	11324202	2
Caelo	Glycerol 85%	11217303	1
Caelo	Glycerol 85%	11355305	1
Caelo	Glycerol 85%	11L06-N02	4
Euro OTC	Glycerol 85%	1201038-01	9
Euro OTC	Glycerol 85%	120103801	1
Gehe	Glycerol 85%	1201038-01	1
Caelo	Glycerol 85%	12069601	2
Caelo	Glycerol 85%	12069604	1
Anzag 17.09.2012	Glycerol 85%	12069603	1
Gehe 08.05.2013	Glycerol 85%	1210047-01	1
Fagron	Glycerol 85%	1209C001	1
Caelo	Glycerol 85%	12069605	1
Noweda	Glycerol 85%	12108201	2
Phönix	Glycerol 85%	12108201	2
Anzag/	Glycerol 85%	12108201	1
Caelo	Glycerol 85%	12108206	1
Phönix 15.12.12	Glycerol 85%	12216601	1
Caelo	Glycerol 85%	12216602	7
Caelo	Glycerol 85%	12216603	17

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	12216604	3
Caelo	Glycerol 85%	12216606	13
GEHE	Glycerol 85%	12216606	2
Gehe 25.04.2013	Glycerol 85%	12216606	1
Caelo	Glycerol 85%	12273604	25
Geo-Chem/sanacorp	Glycerol 85%	12216606	1
Sanacorp	Glycerol 85%	12273604	3
Phönix	Glycerol 85%	12273604	2
Phoenix	Glycerol 85%	12273604	2
A.H. 18.06.13 ek: 6,78EUR	Glycerol 85%	12273604	1
Caelo	Glycerol 85%	12278507	3
Caelo	Glycerol 85%	12278508	4
Fagron	Glycerol 85%	12G11-N02	1
Fagron	Glycerol 85%	12H17-N01	2
Fagron	Glycerol 85%	12J04-N02	1
Caelo	Glycerol 85%	12J04-N09	2
Fagron	Glycerol 85%	12J05-N01	1
Caelo	Glycerol 85%	13039202	6
Caelo	Glycerol 85%	13039203	23
Caelo	Glycerol 85%	13039204	27
Caelo	Glycerol 85%	13060001	5
Caelo	Glycerol 85%	13060006	5
Phönix 14.01.2014	Glycerol 85%	13060001	1
Caelo	Glycerol 85%	13060007	5
Caelo	Glycerol 85%	1309202	1
Fagron	Glycerol 85%	13060006	1
PHOENIX am 14.01.2014	Glycerol 85%	13060007	1
Euro OTC	Glycerol 85%	1306016-02	8
Caelo	Glycerol 85%	13144605	20
Caelo	Glycerol 85%	13144607	24
Gehe	Glycerol 85%	13144607	2
Caelo	Glycerol 85%	13248204	2
Caelo	Glycerol 85%	1327603	1
Fagron	Glycerol 85%	13C28-N02	9
Fagron	Glycerol 85%	13C-N02	1
Fagron	Glycerol 85%	13D11-N05	1
Fagron	Glycerol 85%	13E24-N11	3
Fagron	Glycerol 85%	13J29-B01	1
Caelo	Glycerol 85%	140113PML02	1
Phoenix	Glycerol 85%	1405021-01	1
Laborchemie APOLDA	Glycerol 85%	14149410	1
Caelo	Glycerol 85%	14101304	1
Phönix/	Glycerol 85%	14149410	1
Caelo	Glycerol 85%	16121301	1
Sanacorp	Glycerol 85%		2
Caelo	Glycerol 85%	23081201	1
Caelo	Glycerol 85%	2/AKD52	2
Caelo	Glycerol 85%	2413Q01589	1
Caelo	Glycerol 85%	24021402	1
Caesar/Noweda	Glycerol 85%	2611I-01589	1
Caelo	Glycerol 85%	2405211589	1
Bombastus	Glycerol 85%	276850	2
Bombastus	Glycerol 85%	275075	3
Bombastus	Glycerol 85%	278655	1
Gehe	Glycerol 85%	282731	1
Bombastus	Glycerol 85%	282731	6
Kehr / Lamotte	Glycerol 85%	290114PML04	1
Hedinger/Gehe	Glycerol 85%	3112Q-01589	1
Caelo	Glycerol 85%	4711A-01589	1
Otto Fischar GmbH	Glycerol 85%	5061101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	5x12273604	1
Caelo	Glycerol 85%	5x13039203	1
Kehr / Fischar	Glycerol 85%	5061071	1
Gehe	Glycerol 85%	5062051	1
Fischar/E	Glycerol 85%	5062011	1
Fischar	Glycerol 85%	5062111	5
Otto Fischar GmbH	Glycerol 85%	5062052	4
Caelo	Glycerol 85%	5062052	1
Otto Fischar/Gehe	Glycerol 85%	5063012	1
Phönix	Glycerol 85%	5062111	1
Henry Lamotte GmbH Bremen	Glycerol 85%	5142501	1
Henry Lamotte	Glycerol 85%	5226901	2
Lamotte	Glycerol 85%	5360201	5
Lamotte	Glycerol 85%	5424901	4
Lamotte	Glycerol 85%	5424101	2
Lamotte	Glycerol 85%	5512001	2
Lamotte	Glycerol 85%	5487801	2
Lamotte/Krieger/20.12.13	Glycerol 85%	5538401	4
Caelo	Glycerol 85%		5
HEDINGER	Glycerol 85%	40467	1
Hedinger/Gehe	Glycerol 85%	40483	1
Hedinger	Glycerol 85%	40486	2
HEDINGER	Glycerol 85%	40475	1
hedinger	Glycerol 85%	40494	1
HEDINGER	Glycerol 85%	40497	1
Hedinger/ Sanacorp	Glycerol 85%	40499	1
Hedinger/ Sanacorp	Glycerol 85%	40501	1
HEDINGER	Glycerol 85%	40501	1
Hedinger/Sanacorp	Glycerol 85%	40501	1
Hedinger	Glycerol 85%	40578	1
Caelo	Glycerol 85%	11324208	1
heidringer	Glycerol 85%	40559	1
Caesar & Loretz GmbH	Glycerol 85%	12108201	1
Phoenix	Glycerol 85%	13039203	2
Euro OTC	Glycerol 85%	1306016-01	3
Fagron	Glycerol 85%	1312B03285613	1
Herbeta / Noweda	Glycerol 85%		2
Bombastus	Glycerol 85%	279852	4
Fischar	Glycerol 85%	5063011	3
Spinnrad	Glycerol 85%	713048084	1
Hedinger/Phönix	Glycerol 85%	40531-1	1
Caelo	Glycerol 85%	13136115	1
Hedinger/Phönix	Glycerol 85%	40481	1
Caelo	Glycerol 85%	13276208	5
Anzag,7,17 Euro,16.1014	Glycerol 85%	13356405	1
Caelo	Glycerol 85%	5064011	3
Fischar GmbH/AHD	Glycerol 85%	5064091	1
Otto Fischar	Glycerol 85%	5064021	1
Hedinger Phönix	Glycerol 85%	40507	1
Herbeta Berlin	Glycerol 85%		1
Caelo	Glycerol 85%	P136/14	1
Caelo	Glycerol 85%	13276201	1
Caelo	Glycerol 85%	40480	2
Hedinger Phönix	Glycerol 85%	40495	1
Hedinger Jenne	Glycerol 85%	40501	1
Hedinger/Phönix	Glycerol 85%	40483	2
Hedinger/Sanacorp	Glycerol 85%	40499-3354319071	1
Caelo	Glycerol 85%	1309204	1
Sanacorp	Glycerol 85%	278655	1
Sanacorp	Glycerol 85%	280855	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Glycerol 85%	40475	2
Lamotte/Phönix	Glycerol 85%	124654	2
Henry Lamotte	Glycerol 85%	124654	1
Henry Lamotte Oils GmbH	Glycerol 85%	125545	1
Alliance Health 27.11.2...	Glycerol 85%		1
Bombastus	Glycerol 85%	40215PJK01	1
Sanacorp/ Hedinger	Glycerol 85%	40502	1
Hedinger	Glycerol 85%	40502	1
Hedinger / Alliance	Glycerol 85%	40506	2
Hedinger/ Noweda	Glycerol 85%	40515	1
Hedinger/Gehe	Glycerol 85%	40519	1
Hedinger	Glycerol 85%	40522	1
Aug. Hedinger GmbH & Co	Glycerol 85%	40519	1
Hedinger	Glycerol 85%	40528	1
Caelo	Glycerol 85%	40524	1
Hedinger/Gehe	Glycerol 85%	40524	1
Caelo	Glycerol 85%	40529	1
hedinger/Gehe	Glycerol 85%	40529	1
Hedinger/Noweda	Glycerol 85%	40522	1
Hedinger	Glycerol 85%	40529	1
Hedinger/Gehe	Glycerol 85%	4052B	1
Hedinger/Noweda	Glycerol 85%	40531	1
Hedinger/ Sanacorp	Glycerol 85%	40531-1	1
Hedinger/ Sanacorp	Glycerol 85%	40538	1
Hedinger	Glycerol 85%	40531-1	2
Hedinger	Glycerol 85%	40540	1
hedinger	Glycerol 85%	40538	1
Caelo	Glycerol 85%	40540	1
Hedinger	Glycerol 85%	40541	1
Hedinger/ Sanacorp	Glycerol 85%	40546	2
Hedinger/Noweda	Glycerol 85%	40542	1
Hedinger	Glycerol 85%	40542	1
HirschApotheke	Glycerol 85%	40541	1
Hedinger/Phoenix	Glycerol 85%	40546	1
Hedinger/Phönix	Glycerol 85%	40546	1
Hedinger/Noweda	Glycerol 85%	40557	1
Hedinger/Sanacorp	Glycerol 85%	40557	1
Hedinger/ Sanacorp	Glycerol 85%	40559	1
Hedinger/Noweda	Glycerol 85%	40559	2
Hedinger/Phoenix	Glycerol 85%	40559	1
Hedinger	Glycerol 85%	40559	2
Hedinger/Sanacorp	Glycerol 85%	40566	1
Caelo	Glycerol 85%	40566	2
Hedinger/ Gehe	Glycerol 85%	40566	1
Hedinger/Gehe	Glycerol 85%	40566	1
Stumpf/Hedinger	Glycerol 85%	40576	1
Hedinger	Glycerol 85%	40567	3
Hedinger/Noweda	Glycerol 85%	40566	1
Hedinger/Sanacorp	Glycerol 85%	40576	1
Hedinger	Glycerol 85%	40579	1
aug. hedinger GmbH & Co...	Glycerol 85%	40579	1
Hedinger/Noweda	Glycerol 85%	40583	1
Hedinger/Phönix	Glycerol 85%	40583	1
Hedinger	Glycerol 85%	40583	2
Hedinger	Glycerol 85%	40584	2
Phoenix	Glycerol 85%	40586	1
Aug. Hedinger	Glycerol 85%	40586	1
phönix	Glycerol 85%	40586	1
Hedinger/Gehe	Glycerol 85%	40586	1
phönix	Glycerol 85%	40598	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Noweda	Glycerol 85%	40586	1
Hendiger/Sanacorp	Glycerol 85%	40599	1
Caelo	Glycerol 85%	40598	1
Hedinger/Phönix	Glycerol 85%	40598	1
Hedinger/Phönix	Glycerol 85%	40599	2
Hedinger/Noweda	Glycerol 85%	40599	2
Hedinger/VDL	Glycerol 85%	40599	1
Gehe/Hedinger	Glycerol 85%	40607	1
Hedinger/ Sanacorp	Glycerol 85%	40612	1
Hedinger/Phönix	Glycerol 85%	40612	1
phönix	Glycerol 85%	40612	1
Hedinger/Noweda	Glycerol 85%	40618	1
Gehe - Hedinger	Glycerol 85%	40618	1
Noweda / Hedinger	Glycerol 85%	40626	1
Hedinger/ AHD	Glycerol 85%	40620	1
Hedinger/Sanacorp	Glycerol 85%	40628	1
Hedinger/Phönix	Glycerol 85%	40620	1
Au.HedingerGmbH&CoKG/Al...	Glycerol 85%	40626	1
Hedinger/Phönix	Glycerol 85%	40629	1
Hedinger/Phönix	Glycerol 85%	40628	1
Hedinger/Sanacorp	Glycerol 85%	40629	1
Euro OTC	Glycerol 85%	10022015C	1
Euro OTC	Glycerol 85%	10112014D	1
Caelo	Glycerol 85%	1215302	1
Gehe am 05.04.2014	Glycerol 85%	13276205	2
Cealo	Glycerol 85%	13276205	1
Sanacorp	Glycerol 85%	13276205	1
Caelo	Glycerol 85%	13276206	1
Caelo	Glycerol 85%	13356402	1
	Glycerol 85%	13356404	1
Caelo	Glycerol 85%	13396407	1
AHD	Glycerol 85%	13406202	1
Fagron	Glycerol 85%	13J01-B01-285777	1
Caelo	Glycerol 85%	135406204	1
Fagron	Glycerol 85%	13J01-B01	1
Euro OTC	Glycerol 85%	1401009	1
Caelo	Glycerol 85%	1404024-01	1
Caelo	Glycerol 85%	1405021	1
Euro OTC	Glycerol 85%	140502101	2
Gehe	Glycerol 85%	1405021-02	1
Euro OTC	Glycerol 85%	140903501	1
Caelo	Glycerol 85%	14148104	3
Gehe	Glycerol 85%	14149404	1
sanacorp	Glycerol 85%	14149411	1
Caelo	Glycerol 85%	14149404	1
AHD	Glycerol 85%	14306203	1
Euro OTC	Glycerol 85%	14380504	1
Caelo	Glycerol 85%	14373412	1
Gehe,06.05.15	Glycerol 85%	14306205	1
Fagron	Glycerol 85%	14A10-B01-291983	1
Fagron	Glycerol 85%	14H11-B05-298257	2
Fagron	Glycerol 85%	14D25-B02-293384	1
Caleo	Glycerol 85%	15036903	1
Noweda	Glycerol 85%	1506036-01	1
EuRho/Phönix	Glycerol 85%	1509023-01	2
Caelo	Glycerol 85%	1509023-01	1
Phönix Gotha	Glycerol 85%	15146301	1
Caelo	Glycerol 85%	15146303	1
Caerol/Phoenix	Glycerol 85%	15292203	1
Caesar & Loretz GmbH Hi...	Glycerol 85%	15292204	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
WE:23.03.2016 Sanacorp	Glycerol 85%	15292205	1
Caelo	Glycerol 85%	153782019	1
Caelo	Glycerol 85%	15378201	1
Caelo	Glycerol 85%	15379207	2
Caesar & Loretz/Noweda	Glycerol 85%	15379208	1
Caelo	Glycerol 85%	1601059	1
Caelo	Glycerol 85%	160215301	1
Caelo	Glycerol 85%	160637025	1
Caelo	Glycerol 85%	1607037	1
Euro OTC	Glycerol 85%	160703701	2
Bombastus	Glycerol 85%	1607037-01	1
Fischer/ Noweda	Glycerol 85%	1610032-01	1
Caelo	Glycerol 85%	1610032-01	6
Caelo	Glycerol 85%	161003201	1
Euro OTC	Glycerol 85%	1610032-1	1
GEHE	Glycerol 85%	16215301	1
Caelo	Glycerol 85%	16215308	1
Caelo	Glycerol 85%	16218907	1
Caelo	Glycerol 85%	16218906	2
phönix	Glycerol 85%	16218907	1
Caelo	Glycerol 85%	16297809	1
Caelo	Glycerol 85%	16347206	2
Caelo	Glycerol 85%	16347208	1
Caelo	Glycerol 85%	16355005	1
Henry Lamotte	Glycerol 85%	1703003-01	1
Euro OTC	Glycerol 85%	170509G/1607037-01	1
Euro OTC	Glycerol 85%	1708022-01	2
Caelo	Glycerol 85%	171121G/17157702	1
Caesar und Loretz GmbH	Glycerol 85%	17157702	1
Bombastus	Glycerol 85%	260914PAS01	1
apo ident	Glycerol 85%	1801025-01	1
otto fischer	Glycerol 85%	23081604	1
Caelo	Glycerol 85%	24101401	1
Bombastus	Glycerol 85%	2823731	1
Hedinger	Glycerol 85%	25082015B	1
Bombastus	Glycerol 85%	284401	1
Perrigo/Phönix	Glycerol 85%	284401	1
Bombastus	Glycerol 85%	287120	3
eurho	Glycerol 85%	287120	1
Bombastus	Glycerol 85%	289296	2
Bombastus	Glycerol 85%	290715	3
Bombastus	Glycerol 85%	295502	1
Bombastus	Glycerol 85%	296480	3
Bombastus	Glycerol 85%	298594	1
Bombastus	Glycerol 85%	301816	4
Bombastus	Glycerol 85%	301817	4
Bombastus	Glycerol 85%	302749	3
Euro OTC	Glycerol 85%	302749	2
Bombastus	Glycerol 85%	302750	2
Bombastus	Glycerol 85%	302751	1
Bombastus	Glycerol 85%	305201	4
Bombastus	Glycerol 85%	306830	3
Phönix	Glycerol 85%	40599	1
Bombastus	Glycerol 85%	310593	1
Fischar	Glycerol 85%	5062012	4
Fischar / Phönix	Glycerol 85%	5064011	1
Caelo	Glycerol 85%	5064021	1
Otto Fischer	Glycerol 85%	5064011	1
Alliance Healthcare/Fis...	Glycerol 85%	5063011	1
Otto Fischar GmbH	Glycerol 85%	5064021	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischar	Glycerol 85%	506407	1
Fischar/Phönix	Glycerol 85%	5064071	1
Fischar	Glycerol 85%	5065012	1
Allianc. Healthc./Fischar	Glycerol 85%	5064091	1
fischar/Noweda	Glycerol 85%	5064091	1
Noweda/ Fischar	Glycerol 85%	5065011	1
fishbar	Glycerol 85%	5064091	1
Otto Fischar	Glycerol 85%	5065012	1
Otto Fischar GmbH & Co.KG	Glycerol 85%	5065031	1
Otto Fischar GmbH & Co...	Glycerol 85%	5065012	1
Fischar/Phönix	Glycerol 85%	5065031	1
Fischar	Glycerol 85%	5065061	4
fischar	Glycerol 85%	5065061	2
Caelo	Glycerol 85%	5065061	1
Noweda/ Fischar	Glycerol 85%	5065061	2
Caelo	Glycerol 85%	5065081	1
Fischar/Gehe	Glycerol 85%	5065081	1
Otto Fischer	Glycerol 85%	5065081	1
Fischer	Glycerol 85%	5065101	1
Fischar/ Phönix	Glycerol 85%	5065101	2
Fischar	Glycerol 85%	5065081	1
Gehe/fischar	Glycerol 85%	5065101	1
Fischar	Glycerol 85%	5065101	1
Bombastus	Glycerol 85%	5065101	1
Fischar GmbH/Krieger	Glycerol 85%	5065111	1
otto Fischar	Glycerol 85%	5065111	1
Fischar / Alliance Hea...	Glycerol 85%	5065111	1
Fischar / Noweda	Glycerol 85%	50655081	1
Fischar/ Noweda	Glycerol 85%	5066011	1
Otto Fischar GmbH/ Noweda	Glycerol 85%	5066011	4
Fischar	Glycerol 85%	5066011	1
Otto Fischar GmbH/Gehe	Glycerol 85%	5066011	1
Fischar/Phönix	Glycerol 85%	5066011	1
Fischer	Glycerol 85%	5066021	1
Fischar/ Now	Glycerol 85%	5066031	1
Otto Fischar GmbH & Co.KG	Glycerol 85%	5066021	1
Fischar / Phönix	Glycerol 85%	5066031	1
Fischar GmbH/Krieger	Glycerol 85%	5066021	1
Fischar/Gehe	Glycerol 85%	5066031	1
Fischar	Glycerol 85%	5066031	1
Alliance Healthcare	Glycerol 85%	5066061	1
Otto Fischer GmbH	Glycerol 85%	5066031	1
Otto Fischar GmbH & Co.KG	Glycerol 85%	5066061	1
Otto Fischar	Glycerol 85%	5066031	1
Otto Fischer/Krieger	Glycerol 85%	5066061	1
Otto Fischar	Glycerol 85%	5066061	1
Noweda	Glycerol 85%	5066071	1
Fischar	Glycerol 85%	5066071	1
Fischar/Noweda	Glycerol 85%	5066071	1
fischar/Phönix	Glycerol 85%	5066071	1
Noweda/Otto Fischar GmbH	Glycerol 85%	5066071	1
Fischar/ Now	Glycerol 85%	5066071	4
Fischer	Glycerol 85%	5066071	1
Gehe/Fischar	Glycerol 85%	5066071	1
Fischar/ Phönix	Glycerol 85%	5066081	1
Otto Fischar GmbH & Co.KG	Glycerol 85%	5066111	3
Fischar GmbH/Krieger	Glycerol 85%	5066111	1
Phönix	Glycerol 85%	5066081	1
Caelo	Glycerol 85%	5066081	1
Fischar/Phönix	Glycerol 85%	5066111	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	5066111	1
hochhaus-Apotheke	Glycerol 85%	5066111	1
Otto Fischar	Glycerol 85%	5066111	1
AEP/ Fischar	Glycerol 85%	5067011	1
Otto Fischer	Glycerol 85%	5067011	1
fischar/Phönix	Glycerol 85%	50666061	1
Noweda	Glycerol 85%	5067061	1
Otto Fischar	Glycerol 85%	5067061	1
Sanacorp/Fischar	Glycerol 85%	5067061	1
Otto Fischar GmbH & Co.KG	Glycerol 85%	5067061	1
Fischar/AEP	Glycerol 85%	5067081	1
fischar	Glycerol 85%	5067081	1
Fischar/AHD	Glycerol 85%	5067091	1
Lamotte/Alliance	Glycerol 85%	5603001	1
Henry Lamotte Olis GmbH	Glycerol 85%	5603001	1
Otto Fischar	Glycerol 85%	5068012	1
Henry Lamotte Oils	Glycerol 85%	5538401	6
Henry Lamotte Oils/Phönix	Glycerol 85%	5603001	5
Alliance Health 24.06.2...	Glycerol 85%	5538401	1
Bombastus	Glycerol 85%	6474	1
Sanacorp	Glycerol 85%	72/AKD1	1
Fagron	Glycerol 85%	56-81-5	1
noweda	Glycerol 85%		1
Caelo	Glycerol 85%	73/AKD2	1
Lamotte/AHD	Glycerol 85%	8002208001	1
Euro OTC	Glycerol 85%	7406HH7011	1
Euro OTC	Glycerol 85%	GG61870088	1
Otto	Glycerol 85%		1
Pharmachem GmbH	Glycerol 85%	272305	4
Lamotte/AHD	Glycerol 85%	8001116001	1
Henry Lamotte Phönix	Glycerol 85%	123142	1
Lamotte/Noweda	Glycerol 85%	123464	1
Henry Lamotte	Glycerol 85%	123464	2
Henry Lamotte/Jenne	Glycerol 85%	123464	1
Bartel	Glycerol 85%	40475	2
Heddinge/Sanacorp	Glycerol 85%	40475	2
Hedinger	Glycerol 85%	40478	1
Hedinger	Glycerol 85%	40480	1
Hedinger/Noweda	Glycerol 85%	40480	1
Hedinger/Noweda	Glycerol 85%	40478	2
Hedinger/ Sanacorp	Glycerol 85%	40480	2
hedinger	Glycerol 85%	40483	2
Hedinger/Gehe	Glycerol 85%	40491	1
Hedinger	Glycerol 85%	40481	1
Hedinger/Ebert	Glycerol 85%	40483	1
Gehe	Glycerol 85%	40491	2
Phönix	Glycerol 85%	40481	1
Hedinger	Glycerol 85%	40499	4
Hedinger/Gehe	Glycerol 85%	40497	2
Phönix	Glycerol 85%	40501	1
Noweda/Hedinger	Glycerol 85%	40499	1
hedinger/Gehe	Glycerol 85%	40501	1
Hedinger/Noweda	Glycerol 85%	40502	1
Großhandel	Glycerol 85%	40507	1
Phönix	Glycerol 85%	40507	1
Caelo	Glycerol 85%	40515	1
Hedinger/Gehe	Glycerol 85%	4091	1
Caelo	Glycerol 85%	70913Gly	1
Gehe	Glycerol 85%	1210047-01	1
Euro OTC	Glycerol 85%	1210047-01	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycerol 85%	1213E-01589	1
Gehe 10.04.2013	Glycerol 85%	1227360	1
EuRho	Glycerol 85%	1301025-01	4
Caelo	Glycerol 85%	12278502	1
Caelo	Glycerol 85%	12278501	1
Euro OTC	Glycerol 85%	1301025-01	8
Euro OTC	Glycerol 85%	1301025-02	2
Gehe	Glycerol 85%	1301025-01	1
Caelo	Glycerol 85%	13060002	2
Gehe	Glycerol 85%	1306016-02	1
Caelo	Glycerol 85%	1309203	1
Caesar & Loretz GmbH Hi...	Glycerol 85%	13144607	2
Fagron	Glycerol 85%	13c28-n02	2
Caelo	Glycerol 85%	1529202	1
Caelo	Glycerol 85%	2013A-01589	1
Fagron	Glycerol 85%	250214PJK01	1
Bombastus	Glycerol 85%	276852	4
Gehe	Glycerol 85%	276852	1
Bombastus	Glycerol 85%	2812731	2
Otto Fischar GmbH	Glycerol 85%	5062012	1
Bombastus	Glycerol 85%	280855	1
fischer, Gehe	Glycerol 85%	5062051	1
Fischar/V.d.L.	Glycerol 85%	5062111	1
Fischar/Anzag	Glycerol 85%	5062111	1
Otto Fischer GmbH	Glycerol 85%	5062121	3
Otto Fischar/Phoenix	Glycerol 85%	5063011	1
Gehe	Glycerol 85%	5062121	1
Fischar	Glycerol 85%	5062121	1
Gehe	Glycerol 85%	5063012	1
fischar gmbh	Glycerol 85%	5063012	2
Fischar/Anzag	Glycerol 85%	5063012	1
Lamotte/Noweda	Glycerol 85%	5538401	1
LAMOTTE/NOWEDA	Glycerol 85%	5538401	4
hedinger/anzag	Glycerol 85%	5538401	1
Bombastus	Glycerol 85%	269070	2
Bombastus	Glycerol 85%	266037	1
Gehe	Glycerol 85%	test1	1
Bombastus	Glycerol 85%	271511	2
Gehe	Glycerol 85%	272305	1
Caelo	Glycerol 85%	213E-01589	1
Caelo	Glycerol 85%	3121202	1
Hedinger	Glycerol 85%	40447	1
Hedinger	Glycerol 85%	40452	1
Hedinger / Sanacorp	Glycerol 85%	40451	1
Hedinger/Holdermann 18...	Glycerol 85%	40452	1
Noweda	Glycerol 85%	40464	1
Gehe	Glycerol 85%	40455	1
Phönix 06.10.2012	Glycerol 85%	40467	1
Anzag 20.12.2012	Glycerol 85%	40467	2
Noweda	Glycerol 85%	40473	1
Phönix/Hedinger	Glycerol 85%	40471	1
Sanacorp/ Heddinge	Glycerol 85%	4071	2
Gehe	Glycerol 85%	1001017-01	1
Euro OTC	Glycerol 85%	1108032-01	6
Euro OTC	Glycerol 85%	1108032	1
Caelo	Glycerol 85%	112117302	1
Caelo	Glycerol 85%	1121303	1
Caelo	Glycerol 85%	11324204	1
Caelo	Glycerol 85%	12069603	2
Anzag	Glycerol 85%	120103801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Glycerol 85%	1203C003	1
Caelo	Glycerol 85%	12108202	1
Fagron	Glycerol 85%	12H16-N04	3
Phönix	Glycerol 85%	13039203	1
Euro OTC	Glycerol 85%	12o1038-01	1
Alliance Healthcare	Glycerol 85%	13144605	1
Caelo	Glycerol 85%	19021301	1
Phönix Gotha	Glycerol 85%	275075	1
Henry Lamotte	Glycerol 85%		2
Caelo	Glycerol 85%	20211212-2	1
Noweda	Glycerol 85%	278655	1
Caelo	Glycerol 85%	20061210	1
Bombastus	Glycerol 85%	275076	1
Gehe/Fischer	Glycerol 85%	5060061	1
Fischer otto,Saarbrücken	Glycerol 85%	5060091	1
Hedinger Phönix	Glycerol 85%	40475	1
Otto Fischer GmbH Saarb...	Glycerol 85%	5061011	1
Otto Fischar	Glycerol 85%	5062011	2
Fischar	Glycerol 85%	5061101	1
Fischar	Glycerol 85%	5062051	3
Fischar/Phönix	Glycerol 85%	5062051	1
Otto Fischar	Glycerol 85%	5062051	1
Lamotte	Glycerol 85%	5142501	1
Caelo	Glycerol 85%	5062111	1
Gehe	Glycerol 85%	5226901	1
Noweda	Glycerol 85%	5226901	1
Phönix/Lamotte	Glycerol 85%	5226901	5
Bombastus	Glycerol 85%	266584	1
Fagron	Glycerol 85%	Ch12L04-N09	1
Caesar & Loretz GmbH	Glycerol 85%	12216603	1
Otto Fischar GmbH	Glycerol 85%	5062051	1
HEDINGER	Glycerol 85%	40491	1
Hedinger/Sanacorp	Glycerol 85%	40507	1
Fagron	Glycerol 85%	12L04-N09	1
Sanacorp	Glycerol 85%	13039204	1
Anzag	Glycerol 85%	1306016-02	1
Caesar & Loretz GmbH	Glycerol 85%	13144606	1
Ahd/PharmachemPößneck	Glycerol 85%	13248201	1
Fagron	Glycerol 85%	14A10-B01-291982	1
Fagron	Glycerol 85%	14D25-B02-295895	1
Euro OTC	Glycerol 85%	20122013A	1
Otto Fischar	Glycerol 85%	5060121	1
Caelo	Glycerol 85%	5062121	1
Fischar	Glycerol 85%	5063052	1
Fischar	Glycerol 85%	5063051	1
Hedinger Phönix	Glycerol 85%	40548	1
Hedinger Phönix	Glycerol 85%	40542	1
Hedinger Gehe	Glycerol 85%	40519	1
Anzag,7,62EUR; 27.0415	Glycerol 85%	14306204	1
VDL;16.3.16;7,57EUR	Glycerol 85%	15292205	1
Caelo	Glycerol 85%	1601059-01	1
sanacorp, 7,66EUR, 17.09...	Glycerol 85%	16063701	1
Fischar/AHD	Glycerol 85%	5065031	1
Anzag	Glycerol 85%	40471	1
Anzag	Glycerol 85%	40467	1
Fischar/AHD	Glycerol 85%	5065061	1
Gehe	Glycerol 85%	11355301	1
Gehe	Glycerol 85%	276850	2
Gehe	Glycerol 85%	284401	1

- 89 197 spectra from 501 *Apo-Ident* customers from a total of 26 531 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Glycerol 85%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Glycerol 85%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	430	0	35 114
Type B	0	870	0	92 915
Type C	0	2125	30	89 197

The substance/substance group *Glycerol 85%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.6047 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 99.3103 %)
Type C	100.0000 % (> 99.6586 %)	98.6079 % (> 98.4687 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30886	31737	2.49	58.79
30980	31956	1.54	60.17
31360	31956	1.66	59.41
31377	31737	1.17	61.97
31588	31588	0.00	59.16
31719	31719	0.00	61.70
31737	31737	0.00	62.49
31956	31956	0.00	61.41
34031	34031	0.00	58.83

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Glycolic acid 70%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30491-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Glycolic acid 70%; Acidum glycolicum 70%

Special notes

When selecting the *Glycolic acid 70%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Glycolic acid 70%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Euro OTC	Glycolic acid 70%	1406006-01	31620	40	20140708
Euro OTC	Glycolic acid 70%	1406006-01	31647	40	20140708
Euro OTC	Glycolic acid 70%	1604050	32281	40	20160510*
Euro OTC	Glycolic acid 70%	1604050-03	34184	40	20180315

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Glycolic acid 70%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 384 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 4 reference samples of the substance/substance group *Glycolic acid 70%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Glycolic acid 70%	1201002-01	30491	40
Euro OTC	Glycolic acid 70%	1406006-01	31620 [†]	20
Euro OTC	Glycolic acid 70%	1406006-01	31647 [†]	20
Euro OTC	Glycolic acid 70%	1604050-01	33225	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 13 spectra from 10 *Apo-Ident* customers from 9 batches from the substance/substance group *Glycolic acid 70%*.
- These include spectra of independent samples from 9 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Euro OTC	Glycolic acid 70%	1401014-01	1
Euro OTC	Glycolic acid 70%	1604050-01	3
Euro OTC	Glycolic acid 70%	160405001	2
Euro OTC	Glycolic acid 70%	16040501-01	1
Phoenix	Glycolic acid 70%	1205037-01	1
Euro OTC	Glycolic acid 70%	1604050-02	1
Euro OTC	Glycolic acid 70%	1109027-01	1
Caelo	Glycolic acid 70%	121002201	2
Euro OTC	Glycolic acid 70%	1210022-01	1

- 91 339 spectra from 501 *Apo-Ident* customers from a total of 27 000 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Glycolic acid 70%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Glycolic acid 70%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	120	0	93 754
Type C	0	13	0	91 339

The substance/substance group *Glycolic acid 70%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6592 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31620	32281	1.54	139.87
31647	32281	3.74	134.11
32281	32281	0.00	140.04
34184	32281	8.20	131.95

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Guaiazulene 25%, water-soluble**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30857-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Guaiazulene 25%, water-soluble

Special notes

When selecting the *Guaiazulene 25%, water-soluble* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Guaiazulene 25%, water-soluble*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Guaiazulene 25%,...	12246703	30857	40	20120831
Caelo	Guaiazulene 25%,...	13334401	31157	40	1411565
Caelo	Guaiazulene 25%,...	14130601	31462	40	20140423

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Guaiazulene 25%, water-soluble*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Guaiazulene 25%, water-soluble*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Guaiazulene 25%, water-soluble	14130601	31462 [†]	20
Caelo	Guaiazulene 25%, water-soluble	161905	32321	40
Caelo	Guaiazulene 25%, water-soluble	173365	34100	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 162 spectra from 83 *Apo-Ident* customers from 43 batches from the substance/substance group *Guaiazulene 25%, water-soluble*.
- These include spectra of independent samples from 39 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Sanacorp	Guaiazulene 25%, water-soluble	4/AKD2	2
Phönix	Guaiazulene 25%, water-soluble	12G31-N03	1
Caelo	Guaiazulene 25%, water-soluble	13044906	4
Caelo	Guaiazulene 25%, water-soluble	15075501	9
Caelo	Guaiazulene 25%, water-soluble	16190502	6
Gehe	Guaiazulene 25%, water-soluble	17021503	1
Caelo	Guaiazulene 25%, water-soluble	13334401	6
Euro OTC	Guaiazulene 25%, water-soluble	13334406	1
Caelo	Guaiazulene 25%, water-soluble	516Q-03150	1
Caelo	Guaiazulene 25%, water-soluble	13334407	4
Caelo	Guaiazulene 25%, water-soluble	14130601	9
Caelo	Guaiazulene 25%, water-soluble	14130605	6
Caelo	Guaiazulene 25%, water-soluble	14252104	6
Caelo	Guaiazulene 25%, water-soluble	14252102	6
	Guaiazulene 25%, water-soluble	14B12-B07-292053	1
Fagron	Guaiazulene 25%, water-soluble	14B12-B07-292053	1
sanacorp	Guaiazulene 25%, water-soluble	14G24-B03	1
Fagron	Guaiazulene 25%, water-soluble	14G24B03	1
Fagron	Guaiazulene 25%, water-soluble	14C24-B03-29675	1
anzag	Guaiazulene 25%, water-soluble	14G24-B03-301884	1
Caelo	Guaiazulene 25%, water-soluble	15075503	9
Fagron	Guaiazulene 25%, water-soluble	152716-0001	1
Caelo	Guaiazulene 25%, water-soluble	1515M-03150	1
Caelo	Guaiazulene 25%, water-soluble	15402302	7
Caelo	Guaiazulene 25%, water-soluble	154022302	1
Caelo	Guaiazulene 25%, water-soluble	15402306	6
Caesar & Loretz GmbH Hi...	Guaiazulene 25%, water-soluble	15402306	1
Caelo	Guaiazulene 25%, water-soluble	16190504	2
Caelo	Guaiazulene 25%, water-soluble	17021503	2
Phoenix	Guaiazulene 25%, water-soluble	17021503	1
Fagron	Guaiazulene 25%, water-soluble	17021503	1
gehe	Guaiazulene 25%, water-soluble	17021503	1
Caelo	Guaiazulene 25%, water-soluble	17021508	1
Caelo	Guaiazulene 25%, water-soluble	17192601	3
Caelo	Guaiazulene 25%, water-soluble	17041405	1
Caelo	Guaiazulene 25%, water-soluble	11303702	2
Caelo	Guaiazulene 25%, water-soluble	13334406	3
Caelo	Guaiazulene 25%, water-soluble	13044901	9
Caelo	Guaiazulene 25%, water-soluble	1514U-03150	1
Caelo	Guaiazulene 25%, water-soluble	11303703	1
Fagron	Guaiazulene 25%, water-soluble	11H30-N03	1
Caelo	Guaiazulene 25%, water-soluble	12019802	2
Post-Apotheke	Guaiazulene 25%, water-soluble	12246701	1
Caelo	Guaiazulene 25%, water-soluble	12158701	2
Caelo	Guaiazulene 25%, water-soluble	12246701	3
Sanacorp	Guaiazulene 25%, water-soluble	12246701	2
Caelo	Guaiazulene 25%, water-soluble	12246703	3
Caelo	Guaiazulene 25%, water-soluble	12246704	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Guaiazulene 25%, water-soluble	12G31-N03	1
Noweda	Guaiazulene 25%, water-soluble	12G31-N03	1
Alliance Healthcare	Guaiazulene 25%, water-soluble	13044901	1
Noweda	Guaiazulene 25%, water-soluble	13044901	1
Caelo	Guaiazulene 25%, water-soluble	13044904	8
Großhandel	Guaiazulene 25%, water-soluble	13334401	1
Gehe	Guaiazulene 25%, water-soluble	13334406	1
Fagron	Guaiazulene 25%, water-soluble	13g31-n03	1
Fagron	Guaiazulene 25%, water-soluble	513936	1
Caelo	Guaiazulene 25%, water-soluble		2
Caelo	Guaiazulene 25%, water-soluble	12G31-N03	1
Gehe	Guaiazulene 25%, water-soluble	13334401	1

- 91 190 spectra from 501 *Apo-Ident* customers from a total of 26 968 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Guaiazulene 25%, water-soluble* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Guaiazulene 25%, water-soluble* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	100	0	93 774
Type C	0	150	12	91 190

The substance/substance group *Guaiazulene 25%, water-soluble* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6586 %)	92.5926 % (> 90.7407 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30857	31157	2.72	58.92
31157	31157	0.00	58.54
31462	31157	4.05	61.25

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Hard fat 649**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30770-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Hard fat 649; Adeps solidus 649

Special notes

When selecting the *Hard fat 649* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Hard fat 649*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Hard fat 649	12I10-N12	30770	40	1407601
Fagron	Hard fat 649	15E13-B01-307859	31699	40	20150527

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Hard fat 649*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Hard fat 649*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Hard fat 649	13C19-N03	31118	40
Fagron	Hard fat 649	15E13-B01-307859	31699 [†]	20
Fagron	Hard fat 649	16H30-B03-325618	33249	40
Fagron	Hard fat 649	18H13-B02-352761	34273	40

- 93 734 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 10 spectra from 7 *Apo-Ident* customers from 6 batches from the substance/substance group *Hard fat 649*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Cremer Oleo GmbH	Hard fat 649	150430	1
Cremer Oleo/Anzag	Hard fat 649	150430	1
Fagron	Hard fat 649	15E13-B01-315602	1
Cremer Oleo/Anzag	Hard fat 649	409134	1
Fagron	Hard fat 649	503154	1
Cremer	Hard fat 649	605218	1
Oleochemicals	Hard fat 649	605218	1
Apotheke Leutenberg/AHD	Hard fat 649	605218	1
Apotheke Leutenberg	Hard fat 649	605218	1
Fagron	Hard fat 649	105945	1

- 91 342 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Hard fat 649* can be clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Hard fat 649* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	139	1	93 734
Type C	0	10	0	91 342

The substance/substance group *Hard fat 649* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	99.2857 % (> 97.1429 %)
Type C	100.0000 % (> 99.6594 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to

exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30770	30770	0.00	25.97
31699	30770	10.81	21.22

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Honeys**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30761-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Honeys; Honey; Mel; Mel rosatum; Rose honey

Special notes

When selecting the *Honeys* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Honeys*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Honey	12132005	30761	40	1402562
Caelo	Honey	13269412	31341	40	20130913
Caelo	Rose honey	15406006	31957	40	20160120
Caelo	Rose honey	16104802	32988	40	1708304
Caelo	Rose honey	171208	33333	40	20170619*
Caelo	Honey	171332	33719	40	20170809*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 240 spectra of 6 reference samples from the substance/substance group *Honeys*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 6 different batches.
- 35 304 spectra from a total of 834 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 340 spectra from 9 reference samples of the substance/substance group *Honeys*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Rose honey	12105302	30548	40
Caelo	Honey	11187911	30610	40
Caelo	Rose honey	15406006	31957 [†]	20
Caelo	Honey	150791	32214	40
Caelo	Honey	162659	32698	40
Caelo	Honey	171325	33718	40
Caelo	Honey	171325	33720	40
Caelo	Rose honey	171208	33721	40
Caelo	Rose honey	181215	34050	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 534 spectra from a total of 2159 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 61 spectra from 25 *Apo-Ident* customers from 37 batches from the substance/substance group *Honeys*.
- These include spectra of independent samples from 35 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Honey	13269413	1
Caelo	Honey	13269412	3
Caesar & Loretz GmbH, 4...	Honey	12382405	4
Caesar & Loretz GmbH, 4...	Honey	12132011	2
Caelo	Honey	16265901	1
Caelo	Honey	12382406	1
Euro OTC	Honey	13269401	3
Caesar & Loretz GmbH, 4...	Honey	13269404	1
Caelo	Honey	12132007	1
Caelo	Honey	12132012	1
Caesar & Loretz GmbH, 4...	Honey	12382410	2
Noweda/Apomix	Honey	13269401	1
Caelo	Honey	13269403	4
Phönix	Honey	13269403	1
Caelo	Honey	15079102	1
Caelo	Honey	13269405	2
Noweda	Honey	15079102	1
Caelo	Honey	15079111	1
Noweda	Honey	16265911	1
Caelo	Honey	15079103	1
Caelo	Honey	15417714	1
Caelo	Honey	17132502	2
Caelo	Rose honey	14125607	1
Caelo	Rose honey	5091401	1
Caelo	Rose honey	12105302	1
Caelo	Rose honey	11258901	1
Caelo	Rose honey	12105305	2
Caelo	Rose honey	12315003	1
Caelo	Rose honey	12315005	2
Caelo	Rose honey	13170902	1
Caelo	Rose honey	13170904	1
Caelo	Rose honey	13170906	1
Caelo	Rose honey	15100801	2
Caelo	Rose honey	14125605	1
Caelo	Rose honey	14125603	1
Caelo	Rose honey	15100803	1
Caelo	Rose honey	15406004	2
Caelo	Rose honey	15406006	3
Ebert & Jacobi 02.08.16	Rose honey	15406006	1
Caelo	Rose honey	17120803	1
Caelo	Rose honey	4811M-02884	1

- 91 291 spectra from 501 *Apo-Ident* customers from a total of 26 972 batches from a further 332 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Honeys* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Honeys* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	240	0	35 304
Type B	0	340	0	93 534
Type C	0	60	1	91 291

The substance/substance group *Honeys* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.2353 %)
Type C	100.0000 % (> 99.6587 %)	98.3607 % (> 93.4426 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30761	30761	0.00	36.93
31341	33719	5.03	20.50
31957	32988	2.78	28.10
32988	32988	0.00	26.96
33333	33333	0.00	27.98
33719	33719	0.00	22.51

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Hydrophilic base emulsion (NRF p.25)
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30709-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Hydrophilic base emulsion (NRF p.25); Emulsio basalis hydrophilica

Special notes

When selecting the *Hydrophilic base emulsion (NRF p.25)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Hydrophilic base emulsion (NRF p.25)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Hydrophilic base...	12173002	30709	30	20120628
Caelo	Hydrophilic base...	13305501	31176	40	1403535
Caelo	Hydrophilic base...	14063502	31340	20	20140307
Caelo	Hydrophilic base...	161155	32291	40	20160415*
Caelo	Hydrophilic base...	170444	33154	40	20170217*
Caelo	Hydrophilic base...	163400	33156	40	20161215*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 210 spectra of 6 reference samples from the substance/substance group *Hydrophilic base emulsion (NRF p.25)*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 6 different batches.
- 35 334 spectra from a total of 834 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 671 spectra from 18 reference samples of the substance/substance group *Hydrophilic base emulsion (NRF p.25)*.
- These include spectra of independent samples from 16 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Hydrophilic base emulsion (N...	11129901	30334	41
Caelo	Hydrophilic base emulsion (N...	12173002	30709 [†]	10
Caelo	Hydrophilic base emulsion (N...	14063502	31340 [†]	20
Caelo	Hydrophilic base emulsion (N...	161156	32290	40
Caelo	Hydrophilic base emulsion (N...	161154	32292	40
Caelo	Hydrophilic base emulsion (N...	170443	33155	40
Caelo	Hydrophilic base emulsion (N...	163398	33157	40
Caelo	Hydrophilic base emulsion (N...	163399	33158	40
Caelo	Hydrophilic base emulsion (N...	171184	33412	40
Caelo	Hydrophilic base emulsion (N...	171247	33413	40
Caelo	Hydrophilic base emulsion (N...	171182	33414	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Hydrophilic base emulsion (N...	171183	33415	40
Caelo	Hydrophilic base emulsion (N...	171702	33573	40
Caelo	Hydrophilic base emulsion (N...	171703	33574	40
Caelo	Hydrophilic base emulsion (N...	173421	33971	40
Caelo	Hydrophilic base emulsion (N...	180316	34133	40
Caelo	Hydrophilic base emulsion (N...	181075	34134	40
Caelo	Hydrophilic base emulsion (N...	181806	34248	40

- 93203 spectra from a total of 2149 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1168 spectra from 213 *Apo-Ident* customers from 179 batches from the substance/substance group *Hydrophilic base emulsion (NRF p.25)*.
- These include spectra of independent samples from 173 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic base emulsion (N...	13394302	12
Caelo	Hydrophilic base emulsion (N...	13394304	7
Caelo	Hydrophilic base emulsion (N...	13397601	45
Caelo	Hydrophilic base emulsion (N...	14063502	13
Sanacorp	Hydrophilic base emulsion (N...	14063601	3
Caelo	Hydrophilic base emulsion (N...	14063601	28
Caelo	Hydrophilic base emulsion (N...	14063702	12
Sanacorp	Hydrophilic base emulsion (N...	14141004	1
Caelo	Hydrophilic base emulsion (N...	14238801	16
Caelo	Hydrophilic base emulsion (N...	14238101	1
Caelo	Hydrophilic base emulsion (N...	14238901	28
Caelo	Hydrophilic base emulsion (N...	14140901	12
Caelo	Hydrophilic base emulsion (N...	14239103	8
Caelo	Hydrophilic base emulsion (N...	14360201	9
Caelo	Hydrophilic base emulsion (N...	14351202	6
Caelo	Hydrophilic base emulsion (N...	15209801	22
Caelo	Hydrophilic base emulsion (N...	1515J-03325	1
Caelo	Hydrophilic base emulsion (N...	15210001	10
Caelo	Hydrophilic base emulsion (N...	14360301	29
Caelo	Hydrophilic base emulsion (N...	15439801	20
Caelo	Hydrophilic base emulsion (N...	15337301	19
Caelo	Hydrophilic base emulsion (N...	15440001	18
Caelo	Hydrophilic base emulsion (N...	15440101	20
Caelo	Hydrophilic base emulsion (N...	16115302	9

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic base emulsion (N...	16249301	9
Caelo	Hydrophilic base emulsion (N...	15440102	3
Caelo	Hydrophilic base emulsion (N...	17342101	4
Caelo	Hydrophilic base emulsion (N...	12128401	5
Caelo	Hydrophilic base emulsion (N...	16115303	8
Caelo	Hydrophilic base emulsion (N...	12173002	6
Caelo	Hydrophilic base emulsion (N...	23011503	1
Caelo	Hydrophilic base emulsion (N...	13039701	13
Caelo	Hydrophilic base emulsion (N...	13076702	13
Caelo	Hydrophilic base emulsion (N...	12310101	12
Phönix	Hydrophilic base emulsion (N...	13039701	1
Caelo	Hydrophilic base emulsion (N...	13076803	21
Caelo	Hydrophilic base emulsion (N...	13076801	16
Phönix	Hydrophilic base emulsion (N...	13076803	1
Caelo	Hydrophilic base emulsion (N...	13107703	17
Caelo	Hydrophilic base emulsion (N...	13196001	13
Caelo	Hydrophilic base emulsion (N...	13305401	25
Caelo	Hydrophilic base emulsion (N...	13195801	11
Caelo	Hydrophilic base emulsion (N...	13305601	11
Caelo	Hydrophilic base emulsion (N...	13305501	8
aliance healthcare 5.12...	Hydrophilic base emulsion (N...	13305501	1
Henry Lamotte Oils GmbH...	Hydrophilic base emulsion (N...	13397601	2
Caesar&lorenz,AHC	Hydrophilic base emulsion (N...	14063601	1
Caelo	Hydrophilic base emulsion (N...	14360403	5
Caelo	Hydrophilic base emulsion (N...	16115602	29
Caelo	Hydrophilic base emulsion (N...	30011504	1
Caelo	Hydrophilic base emulsion (N...	28111203	1
Caelo	Hydrophilic base emulsion (N...	30041304	1
Caelo	Hydrophilic base emulsion (N...	914E-03325	1
Caesar & Lorenz	Hydrophilic base emulsion (N...	14238901	1
Caelo	Hydrophilic base emulsion (N...	1514U-03325	1
Caelo	Hydrophilic base emulsion (N...	15337601	5
Caelo	Hydrophilic base emulsion (N...	15209901	10
Caelo	Hydrophilic base emulsion (N...	15337801	20
Caelo/Noweda	Hydrophilic base emulsion (N...	15440001	1
Caelo	Hydrophilic base emulsion (N...	16339901	22
Caelo	Hydrophilic base emulsion (N...	16249201	9
Caelo	Hydrophilic base emulsion (N...	16249302	4
Caelo	Hydrophilic base emulsion (N...	16264101	17
Caelo	Hydrophilic base emulsion (N...	12295002	1
Caelo	Hydrophilic base emulsion (N...	12230203	6
Caelo	Hydrophilic base emulsion (N...	12288103	3
Caelo	Hydrophilic base emulsion (N...	13039702	3
Caelo	Hydrophilic base emulsion (N...	12310102	16
Caelo	Hydrophilic base emulsion (N...	13195901	19
Caesar & Lorenz	Hydrophilic base emulsion (N...	13195901	1
Caelo	Hydrophilic base emulsion (N...	13195902	4
EuRho/Sanac	Hydrophilic base emulsion (N...	13305601	3
Caelo	Hydrophilic base emulsion (N...	1071602	1
Caelo	Hydrophilic base emulsion (N...	3140817	1
Caelo	Hydrophilic base emulsion (N...	2061501	1
Caelo	Hydrophilic base emulsion (N...	11031504	1
Gehe am 6.03.2014	Hydrophilic base emulsion (N...	13195801	1
Caelo	Hydrophilic base emulsion (N...	13305503	6
GEHE 2.04.2014	Hydrophilic base emulsion (N...	13394302	2
Caelo	Hydrophilic base emulsion (N...	13394032	1
Caelo	Hydrophilic base emulsion (N...	13394303	2
caesar & Loretz GmbH	Hydrophilic base emulsion (N...	13394304	2
gehe	Hydrophilic base emulsion (N...	13394402	2
Caelo	Hydrophilic base emulsion (N...	13394402	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar&Loretz GmbH/Anzag	Hydrophilic base emulsion (N...	13397601	1
Caelo	Hydrophilic base emulsion (N...	13730153	1
GEHE 10.03.2014	Hydrophilic base emulsion (N...	14063502	2
AHD	Hydrophilic base emulsion (N...	14063502	2
gehe	Hydrophilic base emulsion (N...	14063503	1
Caelo	Hydrophilic base emulsion (N...	14063503	1
PHOENIX 03.06.2014	Hydrophilic base emulsion (N...	14063701	1
Gehe	Hydrophilic base emulsion (N...	14063702	1
Gehe	Hydrophilic base emulsion (N...	1410901	1
Gehe	Hydrophilic base emulsion (N...	14140901	1
Caelo	Hydrophilic base emulsion (N...	14141004	11
Caelo	Hydrophilic base emulsion (N...	14141003	2
Caelo	Hydrophilic base emulsion (N...	14238701	11
Caelo	Hydrophilic base emulsion (N...	14238702	5
Gehe	Hydrophilic base emulsion (N...	14238702	1
Caelo	Hydrophilic base emulsion (N...	14238802	5
Gehe	Hydrophilic base emulsion (N...	14238802	1
Sanacorp	Hydrophilic base emulsion (N...	14238901	1
Caelo	Hydrophilic base emulsion (N...	14239002	3
Caelo	Hydrophilic base emulsion (N...	1428901	2
Alliance health	Hydrophilic base emulsion (N...	14360201	2
Caelo	Hydrophilic base emulsion (N...	14360202	4
Sanacorp	Hydrophilic base emulsion (N...	14360301	2
Caelo	Hydrophilic base emulsion (N...	14360401	2
Caelo	Hydrophilic base emulsion (N...	1463502	1
Caelo	Hydrophilic base emulsion (N...	15029901	1
Caelo	Hydrophilic base emulsion (N...	15076803	7
Caelo	Hydrophilic base emulsion (N...	15076901	11
Caelo	Hydrophilic base emulsion (N...	15077002	10
Caesar & Loretz GmbH	Hydrophilic base emulsion (N...	15077002	1
Caelo	Hydrophilic base emulsion (N...	15077003	4
Caesar&Loretz GmbH/Anzag	Hydrophilic base emulsion (N...	15077101	2
Caelo	Hydrophilic base emulsion (N...	15077101	21
Cealo	Hydrophilic base emulsion (N...	15077101	3
Caelo	Hydrophilic base emulsion (N...	15081402	1
Noweda	Hydrophilic base emulsion (N...	15209801	1
Hubertus Apotheke	Hydrophilic base emulsion (N...	15209901	1
cealo/noveda	Hydrophilic base emulsion (N...	15209901	1
Caelo	Hydrophilic base emulsion (N...	15209902	6
sanacorp	Hydrophilic base emulsion (N...	15210001	2
Caelo	Hydrophilic base emulsion (N...	15255816	1
Sanacorp	Hydrophilic base emulsion (N...	15337301	2
Caesar&Loretz GmbH	Hydrophilic base emulsion (N...	15337301	1
Bombastus	Hydrophilic base emulsion (N...	15337601	4
Caelo	Hydrophilic base emulsion (N...	15337602	3
Eigen/Defektur	Hydrophilic base emulsion (N...	15337602	1
Caelo	Hydrophilic base emulsion (N...	15337701	6
Caelo	Hydrophilic base emulsion (N...	1537301	1
Caesar & Loretz GmbH Hi...	Hydrophilic base emulsion (N...	15439801	1
Caesar&Loretz GmbH/ Gehe	Hydrophilic base emulsion (N...	15439801	1
Caesar&Loretz GmbH	Hydrophilic base emulsion (N...	15439801	1
Caelo	Hydrophilic base emulsion (N...	15439901	5
Caelo	Hydrophilic base emulsion (N...	15440002	6
Caesar/Noweda	Hydrophilic base emulsion (N...	15440101	1
Caelo	Hydrophilic base emulsion (N...	154401014	1
Caelo	Hydrophilic base emulsion (N...	161115401	1
caesar & loretz	Hydrophilic base emulsion (N...	16115302	1
wepa/Sanaco	Hydrophilic base emulsion (N...	16115303	1
Caelo	Hydrophilic base emulsion (N...	16115401	3
Caesar & Loretz GmbH Hi...	Hydrophilic base emulsion (N...	16115602	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Hydrophilic base emulsion (N...	16115602	1
Ebert 12.09.2016	Hydrophilic base emulsion (N...	16115602	1
Caesar & Loretz GmbH Hi...	Hydrophilic base emulsion (N...	16151602	1
Caelo	Hydrophilic base emulsion (N...	16244101	1
Caesar & Loretz/AHD	Hydrophilic base emulsion (N...	16249201	1
Caelo	Hydrophilic base emulsion (N...	16249202	3
Phoenix	Hydrophilic base emulsion (N...	16249302	1
Noweda	Hydrophilic base emulsion (N...	1626402	1
Caelo	Hydrophilic base emulsion (N...	16264201	4
Allianc. Health. /ca	Hydrophilic base emulsion (N...	16264202	1
Caelo	Hydrophilic base emulsion (N...	16264202	5
Caesar&Loretz GmbH/ All...	Hydrophilic base emulsion (N...	16264202	1
Caelo	Hydrophilic base emulsion (N...	16264203	1
Caelo	Hydrophilic base emulsion (N...	16339801	6
phönix	Hydrophilic base emulsion (N...	16339901	1
Caelo	Hydrophilic base emulsion (N...	16340001	2
Allianz Healthcare	Hydrophilic base emulsion (N...	16340001	1
Caelo	Hydrophilic base emulsion (N...		2
Caelo	Hydrophilic base emulsion (N...	17044401	8
phönix	Hydrophilic base emulsion (N...	17044401	2
Caelo	Hydrophilic base emulsion (N...	17118301	6
Phönix/Phönix	Hydrophilic base emulsion (N...	17118301	1
Sanacorp	Hydrophilic base emulsion (N...	17124701	1
Fagron	Hydrophilic base emulsion (N...	17118401	1
Caelo	Hydrophilic base emulsion (N...	17124701	1
Bayvers	Hydrophilic base emulsion (N...	17124702	2
Caelo	Hydrophilic base emulsion (N...	17124702	7
Caelo	Hydrophilic base emulsion (N...	17170101	1
Caelo	Hydrophilic base emulsion (N...	17170301	5
Caelo	Hydrophilic base emulsion (N...	17170201	1
Caelo	Hydrophilic base emulsion (N...	17226801	6
Caelo	Hydrophilic base emulsion (N...	17226903	2
Caelo	Hydrophilic base emulsion (N...	17342202	1
Caelo	Hydrophilic base emulsion (N...	18031701	1
Caelo	Hydrophilic base emulsion (N...	18081512	1
Caelo	Hydrophilic base emulsion (N...	20061408	1
Caelo	Hydrophilic base emulsion (N...	20061407	1
Caelo	Hydrophilic base emulsion (N...	21031403	1
Caelo	Hydrophilic base emulsion (N...	22042015A	1
Caelo	Hydrophilic base emulsion (N...	23062015D	1
Caelo	Hydrophilic base emulsion (N...	23071507	1
Caelo	Hydrophilic base emulsion (N...	29041406	1
Caelo	Hydrophilic base emulsion (N...	29041502	1
Caelo	Hydrophilic base emulsion (N...	29101402	1
Caelo	Hydrophilic base emulsion (N...	30051501	1
Caelo	Hydrophilic base emulsion (N...	6115401	1
Caelo	Hydrophilic base emulsion (N...	2101304	1
Caelo	Hydrophilic base emulsion (N...	6051304	1
Caelo	Hydrophilic base emulsion (N...	9081301	1
Caelo	Hydrophilic base emulsion (N...	6138389	1
Maros	Hydrophilic base emulsion (N...	5006613	1
Caelo	Hydrophilic base emulsion (N...	10051303	1
Caelo	Hydrophilic base emulsion (N...	10071203	1
Caelo	Hydrophilic base emulsion (N...	10082013A	1
Caelo	Hydrophilic base emulsion (N...	11274003	1
Gehe	Hydrophilic base emulsion (N...	12047602	4
Noweda	Hydrophilic base emulsion (N...	12071901	1
Caelo	Hydrophilic base emulsion (N...	12071901	4
Caelo	Hydrophilic base emulsion (N...	12072002	2
Caesar/Noweda	Hydrophilic base emulsion (N...	12071901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic base emulsion (N...	12121308	1
Gehe	Hydrophilic base emulsion (N...	12072002	1
Gehe	Hydrophilic base emulsion (N...	12128401	1
Sanacorp	Hydrophilic base emulsion (N...	12128401	2
Caelo	Hydrophilic base emulsion (N...	1213003	1
Fiebig ek:4,26EUR 18.09.12	Hydrophilic base emulsion (N...	12173003	1
Caelo	Hydrophilic base emulsion (N...	12173003	1
Caelo	Hydrophilic base emulsion (N...	1218401	1
Anzag 27.11.12 ek:4,33	Hydrophilic base emulsion (N...	12230301	1
Caelo	Hydrophilic base emulsion (N...	12230301	9
Noweda	Hydrophilic base emulsion (N...	12230301	2
Caelo	Hydrophilic base emulsion (N...	12230302	6
Gehe am 07.03.2013	Hydrophilic base emulsion (N...	12310101	1
Caesar & Loretz GmbH, 4...	Hydrophilic base emulsion (N...	12310101	1
Anzag 19.12.12 ek:4,33	Hydrophilic base emulsion (N...	12310102	3
Noweda	Hydrophilic base emulsion (N...	12310102	1
Bombastus	Hydrophilic base emulsion (N...	1288103	1
Phoenix 11.03.13	Hydrophilic base emulsion (N...	13039701	2
A.H. 23.05.13 ek:4,56	Hydrophilic base emulsion (N...	13039701	4
Caelo	Hydrophilic base emulsion (N...	1306801	1
Caesar/Lorenz	Hydrophilic base emulsion (N...	13076702	1
Noweda	Hydrophilic base emulsion (N...	13076801	1
Caelo	Hydrophilic base emulsion (N...	13076802	1
Gehe 29.04.2013	Hydrophilic base emulsion (N...	13076803	2
Celo	Hydrophilic base emulsion (N...	1307681	1
Caelo	Hydrophilic base emulsion (N...	13107702	1
Caesar & Loretz GmbH, 4...	Hydrophilic base emulsion (N...	13107703	1
Caelo	Hydrophilic base emulsion (N...	1310773	2
Gehe	Hydrophilic base emulsion (N...	13195801	1
AHC 15.08.13 EK4,56	Hydrophilic base emulsion (N...	13107703	1
Phönix Gotha	Hydrophilic base emulsion (N...	13195801	1
Caesar & Loretz	Hydrophilic base emulsion (N...	13195801	1
Phönix	Hydrophilic base emulsion (N...	13195801	1
Caesar&Loretz GmbH/Anzag	Hydrophilic base emulsion (N...	13195901	1
Fagron	Hydrophilic base emulsion (N...	13195901	1
Sanacorp	Hydrophilic base emulsion (N...	13195901	1
Euro OTC	Hydrophilic base emulsion (N...	13196001	1
gehe	Hydrophilic base emulsion (N...	13196001	1
Sanacorp	Hydrophilic base emulsion (N...	13305401	5
Fagron	Hydrophilic base emulsion (N...	13305401	5
gehe	Hydrophilic base emulsion (N...	13305401	1
PHOENIX am 27.01.2014	Hydrophilic base emulsion (N...	13305401	1
Eigene Herstellung	Hydrophilic base emulsion (N...	13305601	1
06,15	Hydrophilic base emulsion (N...	13394302	2
Sanacorp	Hydrophilic base emulsion (N...	13397601	2
Anzag	Hydrophilic base emulsion (N...	14063502	3
Phoenix	Hydrophilic base emulsion (N...	14063502	1
Phonix	Hydrophilic base emulsion (N...	14063601	1
Phoenix	Hydrophilic base emulsion (N...	14063701	1
Caelo	Hydrophilic base emulsion (N...	1507702	1
Caelo	Hydrophilic base emulsion (N...	1507002	1
Caelo	Hydrophilic base emulsion (N...	20031306	1
Caelo	Hydrophilic base emulsion (N...	23101202	1
Caelo	Hydrophilic base emulsion (N...	26031305	1
Caelo	Hydrophilic base emulsion (N...	28121206	1
Caelo	Hydrophilic base emulsion (N...	24011403	1
Caelo	Hydrophilic base emulsion (N...	26031303	1
Caelo	Hydrophilic base emulsion (N...	26031302	1
Caesar/Noweda	Hydrophilic base emulsion (N...	5011E-03325	1
Sanacorp	Hydrophilic base emulsion (N...		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Hydrophilic base emulsion (N...	11223501	2
Gehe	Hydrophilic base emulsion (N...	11274003	1
Kehr	Hydrophilic base emulsion (N...	12047601	1
Gehe	Hydrophilic base emulsion (N...	11273902	1
Gehe	Hydrophilic base emulsion (N...	12071902	2
Noweda	Hydrophilic base emulsion (N...	12173003	1
Gehe	Hydrophilic base emulsion (N...	12230203	2
Caelo	Hydrophilic base emulsion (N...	1223030	1
Sanacorp	Hydrophilic base emulsion (N...	12230203	1
Gehe	Hydrophilic base emulsion (N...	13076702	1
Sanacorp	Hydrophilic base emulsion (N...	13039702	1
Phoenix	Hydrophilic base emulsion (N...	13195901	1
Herstellung	Hydrophilic base emulsion (N...	13305501	2
Phoenix	Hydrophilic base emulsion (N...	13397601	1
Caelo	Hydrophilic base emulsion (N...	16071701	1
Caelo	Hydrophilic base emulsion (N...	15061302	1
Caelo	Hydrophilic base emulsion (N...	18111307	1
Caelo	Hydrophilic base emulsion (N...	19091204	1

- 90 184 spectra from 501 *Apo-Ident* customers from a total of 26 837 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Hydrophilic base emulsion (NRF p.25)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Hydrophilic base emulsion (NRF p.25)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	210	0	35 334
Type B	0	671	0	93 203
Type C	0	1126	42	90 183

The substance/substance group *Hydrophilic base emulsion (NRF p.25)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.1429 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 99.1058 %)
Type C	100.0000 % (> 99.6587 %)	96.4041 % (> 96.1473 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30709	33156	3.20	33.70
31176	31176	0.00	32.26
31340	33154	7.93	27.50
32291	32291	0.00	32.00
33154	33154	0.00	34.58
33156	33156	0.00	33.52

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Hydrophilic ointments
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30631-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Hydrophilic ointments; Anionic hydrophilic cream DAB; Anionic hydrophilic cream SR (NRF p.27); Anionic hydrophilic cream with propylene glycol; Anionic hydrophilic cream without preservatives; Base cream, hydrophobic DAC; Cremor anionicus hydrophilicus SR; Cremor basalis hydrophobicus; Cremor nonionicus hydrophilicus SR; Hydrophilic cream, non-ionic SR (NRF p.26); Hydrophilic cream, non-ionic without preservatives SR; Skin care ointment W/L SR; Skin care ointment W/L SR (without preservatives); Transdermal base cream; Unguentum basale hydrophobicum; Unguentum emulsificans aquosum cum propylenglycolum; Unguentum emulsificans aquosum DAB; Unguentum emulsificans aquosum N sine conservante SR; Unguentum emulsificans aquosum N SR (NRF S. 26); Unguentum emulsificans aquosum sine conservante; Unguentum emulsificans aquosum SR (NRF S. 27); Wolff[®] base half-fat cream

Special notes

When selecting the *Hydrophilic ointments* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Hydrophilic ointments*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Dr. Wolff	Wolff® base half. . .	230200	30631	50	20120210
Caelo	Base cream, hydr. . .	12151602	30758	40	20120612
Bombastus	Hydrophilic crea. . .	00273130	30769	40	1405255
Dr. Wolff	Wolff® base half. . .	236420	30837	40	1410330
Caelo	Anionic hydrophi. . .	12221002	30880	50	20120911
apomix	Anionic hydrophi. . .	25A210812	30895	60	20120821
Caelo	Base cream, hydr. . .	12308102	30931	40	20121031
Caelo	Base cream, hydr. . .	13085601	31002	40	20130404
Fagron	Transdermal base. . .	12H17-T02-000281	31008	40	1406052
apomix	Skin care ointme. . .	01B110113	31039	40	1403321
apomix	Skin care ointme. . .	88B111212	31040	40	20121211
Fagron	Transdermal base. . .	12H17-T02-000281	31042	40	20130110
Bombastus	Anionic hydrophi. . .	280574	31054	40	1405103
Dr. Wolff	Wolff® base half. . .	331551	31086	40	20130221
apomix	Hydrophilic crea. . .	21B141013	31188	40	1405386
Caelo	Anionic hydrophi. . .	13343401	31231	40	1403496
Caelo	Anionic hydrophi. . .	13244401	31241	30	20130726
Caelo	Anionic hydrophi. . .	14138901	31440	40	AR-15-FG-008131-01
Caelo	Anionic hydrophi. . .	14063303	31450	40	20140512
Caelo	Anionic hydrophi. . .	14030701	31470	40	20140204
Caelo	Hydrophilic crea. . .	14127002	31477	40	20140512
Caelo	Base cream, hydr. . .	14028101	31478	40	AR-15-FG-009038-01
Dr. Wolff	Wolff® base half. . .	432611	31487	40	20140410
apomix	Hydrophilic crea. . .	14B240614	31500	20	1411613
apomix	Skin care ointme. . .	10A050214	31503	40	1412041
apomix	Anionic hydrophi. . .	11B020914	31537	40	1411614
apomix	Hydrophilic crea. . .	22B061114	31565	40	20141106
Bombastus	Hydrophilic crea. . .	288882	31568	40	20140722
Caelo	Hydrophilic crea. . .	14304902	31591	40	20141013
Fagron	Hydrophilic crea. . .	14E07-M04	31593	40	20140513
apomix	Hydrophilic crea. . .	01B120115	31599	40	20150114
apomix	Skin care ointme. . .	31B020714	31600	40	20140702
Pharmachem	Base cream, hydr. . .	150404	31644	40	20150306
apomix	Anionic hydrophi. . .	05B200415	31661	50	20150422
Caelo	Anionic hydrophi. . .	15048501	31741	40	9128
Caelo	Base cream, hydr. . .	15221202	31799	40	1605111
Caelo	Base cream, hydr. . .	154212	31906	40	20151211*
Caelo	Base cream, hydr. . .	154213	31907	40	20151211*
apomix	Hydrophilic crea. . .	25B081215	31962	40	20151209
Caelo	Hydrophilic crea. . .	153656	31996	40	20151211*
Caelo	Anionic hydrophi. . .	154004	32022	30	20151117*
Caelo	Anionic hydrophi. . .	153971	32039	40	20151113*
Caelo	Anionic hydrophi. . .	153573	32052	40	20151015*
Caelo	Anionic hydrophi. . .	153295	32123	40	20150924*
Caelo	Anionic hydrophi. . .	160673	32132	40	20160310*
Caelo	Anionic hydrophi. . .	160677	32142	40	20160310*
Pharmachem	Base cream, hydr. . .	160825	32303	40	20160517
Fagron	Transdermal base. . .	15A28-T10-018922	32346	40	20150804
Caelo	Anionic hydrophi. . .	161525	32435	40	20160518*
Caelo	Anionic hydrophi. . .	161916	32620	40	20160630*
Caelo	Anionic hydrophi. . .	161917	32621	40	20160630*
Caelo	Anionic hydrophi. . .	162901	32882	40	20161018*
Pharmachem	Skin care ointme. . .	170117	33295	40	20180131
Caelo	Anionic hydrophi. . .	171023	33354	40	20170419*
Caelo	Base cream, hydr. . .	171920	33595	40	20170710*
apomix	Skin care ointme. . .	24B200617	33641	40	20170620
Fagron	Base cream, hydr. . .	17C06-T08-042287	33964	40	20170425

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 2290 spectra of 57 reference samples from the substance/substance group *Hydrophilic ointments*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 56 different batches.
- 33 254 spectra from a total of 784 batches from further 201 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 14 173 spectra from 343 reference samples of the substance/substance group *Hydrophilic ointments*.
- These include spectra of independent samples from 223 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Bombastus	Base cream, hydrophobic DAC	00234743	30219	40
zu erfassen	Anionic hydrophilic cream wit...	zu erfassen	30223	40
Bombastus	Hydrophilic cream, non-ionic...	00258017	30233	40
Bombastus	Anionic hydrophilic cream DAB	00263988	30333	40
Bombastus	Anionic hydrophilic cream SR...	00262477	30337	40
apomix	Anionic hydrophilic cream SR...	24A080811	30349	20
apomix	Anionic hydrophilic cream wit...	15B111011	30350	30
apomix	Hydrophilic cream, non-ionic...	25A140911	30352	30
apomix	Hydrophilic cream, non-ionic...	27B111011	30353	30
apomix	Skin care ointment W/L SR	58B190711	30355	30
Leyh Pharma	Hydrophilic cream, non-ionic...	112311	30368	30
apomix	Hydrophilic cream, non-ionic...	01A030112	30392	40
Bombastus	Anionic hydrophilic cream DAB	00266260	30402	40
Caelo	Anionic hydrophilic cream DAB	11225101	30404	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Hydrophilic cream, non-ionic...	34B191211	30440	40
Leyh Pharma	Hydrophilic cream, non-ionic...	120211	30441	40
apomix	Hydrophilic cream, non-ionic...	00267186	30444	40
Caelo	Anionic hydrophilic cream DAB	11258201	30446	40
Caelo	Anionic hydrophilic cream DAB	11326601	30462	40
Caelo	Base cream, hydrophobic DAC	12025003	30467	60
Fagron	Hydrophilic cream, non-ionic...	09B190312	30483	40
Caelo	Anionic hydrophilic cream wi...	12117005	30572	40
Caelo	Anionic hydrophilic cream wi...	12117002	30577	80
Dr. Wolff	Wolff® base half-fat cream	zu erfassen	30715	40
apomix	Skin care ointment W/L SR	10A050214	31503 [†]	20
apomix	Hydrophilic cream, non-ionic...	14B240614	31516	60
apomix	Anionic hydrophlic cream wit...	11B020914	31537 [†]	20
apomix	Hydrophilic cream, non-ionic...	22B061114	31565 [†]	20
Bombastus	Hydrophilic cream, non-ionic...	288882	31568 [†]	20
Caelo	Hydrophilic cream, non-ionic...	14304902	31591 [†]	20
Fagron	Hydrophilic cream, non-ionic...	14E07-M04	31593 [†]	20
apomix	Hydrophilic cream, non-ionic...	01B120115	31599 [†]	20
apomix	Skin care ointment W/L SR (w...	31B020714	31600 [†]	20
Pharmachem	Base cream, hydrophobic DAC	150404	31644 [†]	20
apomix	Anionic hydrophlic cream wit...	05B200415	31661 [†]	20
Caelo	Anionic hydrophilic cream wi...	15048501	31741 [†]	20
Caelo	Base cream, hydrophobic DAC	15221202	31799 [†]	20
Caelo	Base cream, hydrophobic DAC	154212	31906 [†]	20
Caelo	Base cream, hydrophobic DAC	154213	31907 [†]	20
Caelo	Base cream, hydrophobic DAC	154214	31908	60
apomix	Hydrophilic cream, non-ionic...	25B081215	31962 [†]	20
Caelo	Hydrophilic cream, non-ionic...	160175	31994	60
Caelo	Hydrophilic cream, non-ionic...	153258	31995	60
Caelo	Anionic hydrophilic cream DAB	154004	32022 [†]	15
Caelo	Anionic hydrophilic cream DAB	154005	32023	60
Caelo	Anionic hydrophilic cream DAB	154007	32024	60
Caelo	Anionic hydrophilic cream DAB	154009	32025	60
Caelo	Anionic hydrophilic cream DAB	154012	32026	60
Caelo	Anionic hydrophilic cream DAB	154012	32027	60
Caelo	Anionic hydrophilic cream DAB	154013	32028	60
Caelo	Anionic hydrophilic cream DAB	154014	32029	60
Caelo	Anionic hydrophilic cream DAB	154015	32030	60
Caelo	Anionic hydrophilic cream DAB	154317	32031	60
Caelo	Anionic hydrophilic cream DAB	154317	32032	60
Caelo	Anionic hydrophilic cream DAB	154321	32033	60
Caelo	Anionic hydrophilic cream DAB	154321	32034	60
Caelo	Anionic hydrophilic cream DAB	153529	32035	60
Caelo	Anionic hydrophilic cream DAB	153530	32036	60
Caelo	Anionic hydrophilic cream DAB	153970	32037	60
Caelo	Anionic hydrophilic cream DAB	153971	32038	60
Caelo	Anionic hydrophilic cream DAB	153971	32039 [†]	20
Caelo	Anionic hydrophilic cream DAB	153972	32040	60
Caelo	Anionic hydrophilic cream DAB	153973	32041	60
Caelo	Anionic hydrophilic cream DAB	153575	32042	60
Caelo	Anionic hydrophilic cream DAB	153969	32043	60
Caelo	Anionic hydrophilic cream DAB	153969	32044	60
Caelo	Anionic hydrophilic cream DAB	154281	32045	60
Caelo	Anionic hydrophilic cream DAB	154282	32046	60
Caelo	Anionic hydrophilic cream DAB	154282	32047	60
Caelo	Anionic hydrophilic cream DAB	153531	32048	60
Caelo	Anionic hydrophilic cream DAB	153571	32049	60
Caelo	Anionic hydrophilic cream DAB	153572	32050	60
Caelo	Anionic hydrophilic cream DAB	153574	32051	60

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Anionic hydrophilic cream DAB	153573	32052 [†]	20
Caelo	Anionic hydrophilic cream DAB	153573	32053	60
Caelo	Anionic hydrophilic cream DAB	153575	32054	60
Caelo	Anionic hydrophilic cream DAB	154284	32055	60
Caelo	Anionic hydrophilic cream DAB	154284	32056	60
Caelo	Anionic hydrophilic cream DAB	153194	32057	60
Caelo	Anionic hydrophilic cream DAB	154281	32058	60
Caelo	Anionic hydrophilic cream DAB	154285	32059	60
Caelo	Anionic hydrophilic cream DAB	154285	32060	60
Caelo	Anionic hydrophilic cream DAB	154283	32061	60
Caelo	Anionic hydrophilic cream DAB	154283	32062	60
Caelo	Anionic hydrophilic cream DAB	153968	32063	60
Caelo	Anionic hydrophilic cream SR...	161075	32115	40
Caelo	Anionic hydrophilic cream SR...	161075	32116	40
Caelo	Hydrophilic cream, non-ionic...	161141	32119	40
Caelo	Hydrophilic cream, non-ionic...	161141	32120	40
Caelo	Hydrophilic cream, non-ionic...	161142	32121	40
Caelo	Hydrophilic cream, non-ionic...	161142	32122	40
Caelo	Anionic hydrophilic cream SR...	153296	32124	40
Caelo	Anionic hydrophilic cream SR...	153865	32125	40
Caelo	Anionic hydrophilic cream SR...	153866	32126	40
Caelo	Anionic hydrophilic cream SR...	160196	32127	40
Caelo	Anionic hydrophilic cream SR...	160198	32128	40
Caelo	Anionic hydrophilic cream wi...	160673	32129	40
Caelo	Anionic hydrophilic cream wi...	160673	32130	40
Caelo	Anionic hydrophilic cream wi...	160673	32131	40
Caelo	Anionic hydrophilic cream wi...	160674	32133	40
Caelo	Anionic hydrophilic cream wi...	160674	32134	40
Caelo	Anionic hydrophilic cream wi...	160674	32135	40
Caelo	Anionic hydrophilic cream wi...	160674	32136	40
Caelo	Anionic hydrophilic cream wi...	160675	32137	40
Caelo	Anionic hydrophilic cream wi...	160675	32138	40
Caelo	Anionic hydrophilic cream wi...	160676	32139	40
Caelo	Anionic hydrophilic cream wi...	160676	32140	40
Caelo	Anionic hydrophilic cream wi...	160677	32141	40
Caelo	Anionic hydrophilic cream wi...	160678	32143	40
Caelo	Anionic hydrophilic cream wi...	160679	32144	40
Caelo	Anionic hydrophilic cream DAB	160578	32145	40
Caelo	Anionic hydrophilic cream DAB	161031	32146	40
Caelo	Anionic hydrophilic cream DAB	161031	32147	40
Caelo	Anionic hydrophilic cream DAB	161033	32148	40
Caelo	Base cream, hydrophobic DAC	160748	32149	40
Caelo	Base cream, hydrophobic DAC	160748	32150	40
Caelo	Base cream, hydrophobic DAC	160749	32151	40
Caelo	Base cream, hydrophobic DAC	160749	32152	40
Caelo	Base cream, hydrophobic DAC	160751	32153	40
Caelo	Base cream, hydrophobic DAC	160751	32154	40
Caelo	Base cream, hydrophobic DAC	160750	32155	40
Caelo	Base cream, hydrophobic DAC	160752	32156	40
Caelo	Anionic hydrophilic cream DAB	161510	32218	40
Caelo	Anionic hydrophilic cream DAB	161510	32219	40
Caelo	Anionic hydrophilic cream DAB	161507	32220	40
Caelo	Anionic hydrophilic cream DAB	161507	32221	40
Caelo	Anionic hydrophilic cream DAB	161508	32222	40
Caelo	Anionic hydrophilic cream SR...	161267	32225	40
Caelo	Base cream, hydrophobic DAC	160750	32226	40
Caelo	Base cream, hydrophobic DAC	161654	32370	41
Caelo	Base cream, hydrophobic DAC	161654	32371	41
Caelo	Base cream, hydrophobic DAC	161655	32372	41
Caelo	Base cream, hydrophobic DAC	161655	32373	41

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Base cream, hydrophobic DAC	161656	32374	41
Caelo	Base cream, hydrophobic DAC	161656	32375	40
Caelo	Base cream, hydrophobic DAC	161657	32376	41
Caelo	Base cream, hydrophobic DAC	161657	32377	41
Caelo	Base cream, hydrophobic DAC	161742	32378	41
Caelo	Hydrophilic cream, non-ionic...	161694	32382	40
Caelo	Hydrophilic cream, non-ionic...	161694	32383	40
Caelo	Hydrophilic cream, non-ionic...	161695	32384	40
Caelo	Hydrophilic cream, non-ionic...	161695	32385	40
Caelo	Anionic hydrophilic cream SR...	162155	32386	40
Caelo	Anionic hydrophilic cream SR...	162156	32387	40
Caelo	Anionic hydrophilic cream SR...	162157	32388	40
Caelo	Anionic hydrophilic cream wi...	161918	32422	40
Caelo	Anionic hydrophilic cream wi...	161918	32423	40
Caelo	Anionic hydrophilic cream wi...	161919	32424	40
Caelo	Anionic hydrophilic cream wi...	161919	32425	40
Caelo	Anionic hydrophilic cream DAB	161506	32426	40
Caelo	Anionic hydrophilic cream DAB	161506	32427	40
Caelo	Anionic hydrophilic cream DAB	161506	32428	40
Caelo	Anionic hydrophilic cream DAB	161510	32429	40
Caelo	Anionic hydrophilic cream DAB	161510	32430	40
Caelo	Anionic hydrophilic cream DAB	161798	32431	40
Caelo	Anionic hydrophilic cream DAB	161798	32432	40
Caelo	Anionic hydrophilic cream DAB	161798	32433	40
Caelo	Anionic hydrophilic cream DAB	161798	32434	40
Caelo	Anionic hydrophilic cream DAB	161525	32436	40
Caelo	Anionic hydrophilic cream DAB	161529	32437	40
Caelo	Anionic hydrophilic cream DAB	161529	32438	40
Caelo	Anionic hydrophilic cream DAB	161529	32439	40
Caelo	Anionic hydrophilic cream DAB	161526	32440	40
Caelo	Anionic hydrophilic cream DAB	161526	32441	40
Caelo	Anionic hydrophilic cream DAB	161526	32442	40
Caelo	Anionic hydrophilic cream DAB	161526	32443	40
Caelo	Anionic hydrophilic cream DAB	161527	32444	40
Caelo	Anionic hydrophilic cream DAB	161527	32445	40
Caelo	Anionic hydrophilic cream DAB	161527	32446	40
Caelo	Anionic hydrophilic cream DAB	161528	32447	40
Caelo	Anionic hydrophilic cream DAB	161528	32448	40
Caelo	Anionic hydrophilic cream DAB	161528	32449	40
Caelo	Anionic hydrophilic cream DAB	161560	32450	40
Caelo	Anionic hydrophilic cream DAB	161561	32451	40
Caelo	Anionic hydrophilic cream DAB	161563	32452	40
Caelo	Anionic hydrophilic cream DAB	161564	32453	40
Caelo	Base cream, hydrophobic DAC	161654	32491	40
Caelo	Base cream, hydrophobic DAC	161655	32492	40
Caelo	Base cream, hydrophobic DAC	161742	32508	40
Caelo	Base cream, hydrophobic DAC	161656	32509	40
Caelo	Base cream, hydrophobic DAC	161657	32510	40
Caelo	Anionic hydrophilic cream SR...	162586	32550	40
Caelo	Anionic hydrophilic cream SR...	162587	32551	40
Caelo	Anionic hydrophilic cream SR...	162587	32552	40
Caelo	Anionic hydrophilic cream SR...	162587	32553	40
Caelo	Hydrophilic cream, non-ionic...	161694	32582	40
Caelo	Hydrophilic cream, non-ionic...	161695	32583	40
Caelo	Hydrophilic cream, non-ionic...	162220	32584	40
Caelo	Hydrophilic cream, non-ionic...	162220	32585	40
Caelo	Hydrophilic cream, non-ionic...	162220	32594	40
Caelo	Hydrophilic cream, non-ionic...	162220	32595	40
Caelo	Hydrophilic cream, non-ionic...	162221	32596	40
Caelo	Hydrophilic cream, non-ionic...	162221	32597	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Hydrophilic cream, non-ionic...	162221	32598	40
Caelo	Hydrophilic cream, non-ionic...	162221	32599	40
Caelo	Hydrophilic cream, non-ionic...	162221	32600	40
Caelo	Hydrophilic cream, non-ionic...	162221	32601	40
Caelo	Hydrophilic cream, non-ionic...	162579	32602	40
Caelo	Hydrophilic cream, non-ionic...	162579	32603	40
Caelo	Anionic hydrophilic cream wi...	161915	32619	40
Caelo	Anionic hydrophilic cream wi...	161918	32622	40
Caelo	Anionic hydrophilic cream wi...	161919	32623	40
Caelo	Anionic hydrophilic cream DAB	161561	32642	40
Caelo	Anionic hydrophilic cream DAB	161563	32643	40
Caelo	Anionic hydrophilic cream DAB	161589	32644	40
Caelo	Anionic hydrophilic cream DAB	161796	32645	40
Caelo	Anionic hydrophilic cream DAB	161796	32646	40
Caelo	Anionic hydrophilic cream DAB	161797	32647	40
Caelo	Anionic hydrophilic cream DAB	161800	32648	40
Caelo	Anionic hydrophilic cream DAB	162500	32649	40
Caelo	Anionic hydrophilic cream DAB	162500	32650	40
Caelo	Anionic hydrophilic cream DAB	162501	32651	40
Caelo	Anionic hydrophilic cream DAB	162501	32652	40
Caelo	Anionic hydrophilic cream DAB	162502	32653	40
Caelo	Anionic hydrophilic cream DAB	162589	32654	40
Caelo	Anionic hydrophilic cream DAB	162589	32655	40
Caelo	Anionic hydrophilic cream DAB	162590	32656	40
Caelo	Anionic hydrophilic cream DAB	162590	32657	40
Caelo	Anionic hydrophilic cream DAB	162591	32658	40
Caelo	Anionic hydrophilic cream DAB	162591	32659	40
Caelo	Anionic hydrophilic cream DAB	162591	32660	40
Caelo	Anionic hydrophilic cream DAB	162592	32661	40
Caelo	Anionic hydrophilic cream DAB	162592	32662	40
Caelo	Anionic hydrophilic cream DAB	162593	32663	40
Caelo	Anionic hydrophilic cream DAB	162593	32664	40
Caelo	Anionic hydrophilic cream DAB	162593	32665	40
Caelo	Anionic hydrophilic cream DAB	162593	32666	40
Caelo	Anionic hydrophilic cream DAB	162594	32667	40
Caelo	Anionic hydrophilic cream DAB	162594	32668	40
Caelo	Anionic hydrophilic cream DAB	162596	32669	40
Caelo	Anionic hydrophilic cream DAB	162596	32670	40
Caelo	Anionic hydrophilic cream DAB	162596	32671	40
Caelo	Anionic hydrophilic cream DAB	162597	32672	40
Caelo	Anionic hydrophilic cream DAB	162597	32673	40
Caelo	Anionic hydrophilic cream DAB	162599	32674	40
Caelo	Anionic hydrophilic cream DAB	162599	32675	40
Caelo	Anionic hydrophilic cream DAB	162600	32676	40
Caelo	Anionic hydrophilic cream DAB	162600	32677	40
Caelo	Anionic hydrophilic cream DAB	162645	32678	40
Caelo	Base cream, hydrophobic DAC	162937	32855	40
Caelo	Base cream, hydrophobic DAC	162938	32856	40
Caelo	Base cream, hydrophobic DAC	162939	32857	40
Caelo	Anionic hydrophilic cream DAB	162643	32875	40
Caelo	Anionic hydrophilic cream DAB	162644	32876	40
Caelo	Anionic hydrophilic cream DAB	162646	32877	40
Caelo	Anionic hydrophilic cream DAB	162647	32878	40
Caelo	Anionic hydrophilic cream DAB	162598	32879	40
Caelo	Anionic hydrophilic cream DAB	163172	32880	40
Caelo	Anionic hydrophilic cream DAB	163226	32881	40
Caelo	Anionic hydrophilic cream wi...	162902	32883	40
Caelo	Anionic hydrophilic cream wi...	162903	32884	40
Caelo	Anionic hydrophilic cream wi...	162904	32885	40
Caelo	Anionic hydrophilic cream wi...	162905	32886	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Anionic hydrophilic cream DAB	170451	33122	40
Caelo	Anionic hydrophilic cream DAB	170452	33123	40
Caelo	Anionic hydrophilic cream SR...	170419	33125	40
Caelo	Base cream, hydrophobic DAC	170732	33127	40
Caelo	Anionic hydrophilic cream wi...	170302	33128	40
Caelo	Anionic hydrophilic cream wi...	170303	33129	40
Caelo	Anionic hydrophilic cream wi...	170299	33136	40
Caelo	Anionic hydrophilic cream wi...	170300	33137	40
Caelo	Anionic hydrophilic cream wi...	170301	33138	40
Caelo	Anionic hydrophilic cream SR...	170418	33139	40
Caelo	Anionic hydrophilic cream DAB	163281	33170	40
Caelo	Anionic hydrophilic cream DAB	163282	33171	40
Caelo	Anionic hydrophilic cream DAB	163283	33172	40
Caelo	Anionic hydrophilic cream DAB	163284	33173	40
Caelo	Anionic hydrophilic cream DAB	170453	33174	40
Caelo	Anionic hydrophilic cream DAB	170454	33175	40
Caelo	Hydrophilic cream, non-ionic...	170097	33187	40
apomix	Skin care ointment W/L SR (w...	14B060417	33264	40
apomix	Hydrophilic cream, non-ionic...	07B030417	33283	40
Pharmachem	Skin care ointment W/L SR	180117	33294	40
Caelo	Anionic hydrophilic cream SR...	171258	33338	40
Caelo	Anionic hydrophilic cream SR...	171442	33339	40
Caelo	Anionic hydrophilic cream wi...	171021	33350	40
Caelo	Anionic hydrophilic cream wi...	171024	33351	40
Caelo	Anionic hydrophilic cream wi...	170302	33352	40
Caelo	Anionic hydrophilic cream wi...	171020	33353	40
Caelo	Anionic hydrophilic cream wi...	171022	33355	40
Caelo	Anionic hydrophilic cream wi...	170299	33356	40
Caelo	Base cream, hydrophobic DAC	171504	33373	40
Caelo	Base cream, hydrophobic DAC	171505	33374	40
Caelo	Base cream, hydrophobic DAC	171130	33375	30
Caelo	Base cream, hydrophobic DAC	171129	33376	30
Caelo	Base cream, hydrophobic DAC	171504	33377	30
Caelo	Hydrophilic cream, non-ionic...	171134	33378	40
Caelo	Hydrophilic cream, non-ionic...	171140	33379	40
Fagron	Transdermal base cream	1609813	33479	40
Dr. Wolff	Wolff® base half-fat cream	702500	33484	40
Dr. Wolff	Wolff® base half-fat cream	701451	33485	40
apomix	Anionic hydrophilic cream wit...	05B310517	33516	40
apomix	Anionic hydrophilic cream wit...	05B310517	33517	40
Caelo	Hydrophilic cream, non-ionic...	172126	33551	40
Caelo	Anionic hydrophilic cream SR...	172117	33552	40
Caelo	Anionic hydrophilic cream SR...	172118	33553	40
Caelo	Anionic hydrophilic cream wi...	171565	33557	40
Caelo	Anionic hydrophilic cream wi...	171566	33558	40
Caelo	Anionic hydrophilic cream wi...	171854	33559	40
Caelo	Anionic hydrophilic cream DAB	171421	33581	40
Caelo	Anionic hydrophilic cream DAB	171422	33582	40
Caelo	Anionic hydrophilic cream DAB	171423	33583	40
Caelo	Anionic hydrophilic cream DAB	171424	33584	40
Caelo	Anionic hydrophilic cream DAB	171427	33585	40
Caelo	Anionic hydrophilic cream DAB	171430	33586	40
Caelo	Anionic hydrophilic cream DAB	171431	33587	40
Caelo	Anionic hydrophilic cream DAB	171432	33588	40
Caelo	Anionic hydrophilic cream DAB	171433	33589	40
Caelo	Anionic hydrophilic cream DAB	171434	33590	40
Caelo	Anionic hydrophilic cream DAB	171435	33591	40
Caelo	Anionic hydrophilic cream DAB	171436	33592	30
Caelo	Anionic hydrophilic cream DAB	171437	33593	40
Caelo	Anionic hydrophilic cream DAB	171438	33594	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
apomix	Skin care ointment W/L SR	31A130717	33627	40
apomix	Hydrophilic cream, non-ionic...	16B230817	33643	40
Caelo	Hydrophilic cream, non-ionic...	172610	33754	40
Caelo	Anionic hydrophilic cream wi...	172621	33755	40
Caelo	Anionic hydrophilic cream wi...	172619	33756	40
Caelo	Anionic hydrophilic cream wi...	172622	33757	40
Caelo	Anionic hydrophilic cream SR...	172568	33768	40
Caelo	Anionic hydrophilic cream SR...	172569	33769	40
apomix	Skin care ointment W/L SR (w...	37A021117	33792	40
apomix	Anionic hydrophilic cream wit...	01B080118	33947	40
Caelo	Hydrophilic cream, non-ionic...	173139	34019	40
Caelo	Hydrophilic cream, non-ionic...	172940	34020	40
Caelo	Anionic hydrophilic cream SR...	172944	34021	40
Caelo	Anionic hydrophilic cream SR...	173195	34022	40
Caelo	Base cream, hydrophobic DAC	172552	34023	40
Caelo	Anionic hydrophilic cream wi...	173353	34028	40
Caelo	Anionic hydrophilic cream wi...	180414	34140	40
Caelo	Hydrophilic cream, non-ionic...	180503	34144	40
Caelo	Hydrophilic cream, non-ionic...	181752	34145	40
Caelo	Anionic hydrophilic cream SR...	181023	34152	40
Caelo	Anionic hydrophilic cream SR...	181475	34153	40
Caelo	Anionic hydrophilic cream SR...	180411	34162	40
Caelo	Base cream, hydrophobic DAC	180512	34163	40
Caelo	Base cream, hydrophobic DAC	181739	34164	40
Fagron	Transdermal base cream	1802650/2	34204	40
Caelo	Anionic hydrophilic cream SR...	182181	34260	40
Caelo	Anionic hydrophilic cream wi...	181808	34262	40
Caelo	Anionic hydrophilic cream DAB	182434	34264	40
Caelo	Base cream, hydrophobic DAC	182398	34269	40
Caelo	Hydrophilic cream, non-ionic...	182180	34272	40

- 79 701 spectra from a total of 1922 batches from further 426 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 9850 spectra from 457 *Apo-Ident* customers from 2982 batches from the substance/substance group *Hydrophilic ointments*.
- These include spectra of independent samples from 2933 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
AMH Niemann / E	Anionic hydrophilic cream SR...	1A050116	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Anionic hydrophilic cream SR...	1A020117	1
AMH Niemann/G	Anionic hydrophlic cream wit...	1B070114	1
Apomix	Anionic hydrophilic cream SR...	1A050116	1
apomix/Phönix	Anionic hydrophlic cream wit...	1B160117	1
AMH Niemann / Noweda	Anionic hydrophlic cream wit...	1B070114	1
Apomix	Anionic hydrophlic cream wit...	1B160117	1
Caelo	Anionic hydrophilic cream SR...	2091403	1
Bombastus	Anionic hydrophilic cream SR...	292980	2
AMH Niemann / G	Anionic hydrophilic cream SR...	2A130114	1
PKH GmbH	Anionic hydrophilic cream SR...	2A030117	1
PKH Halle	Anionic hydrophlic cream wit...	2B060317	2
Gehe	Anionic hydrophlic cream wit...	2B100214	1
Gehe	Anionic hydrophilic cream SR...	3A120116	2
apomix/Noweda	Anionic hydrophlic cream wit...	2B250116	1
Apomix/Noweda	Anionic hydrophilic cream SR...	3A120116	2
AMH Niemann / G	Anionic hydrophilic cream SR...	3A200115	1
apomix/Phönix	Anionic hydrophilic cream SR...	3A310117	1
AMH Niemann / K	Anionic hydrophilic cream SR...	3A310117	1
Apomix/AHD	Anionic hydrophlic cream wit...	3B170216	1
phönix	Anionic hydrophilic cream SR...	3a310117	1
AMH Niemann GmbH&Co.KG	Anionic hydrophlic cream wit...	3B020315	1
Caelo	Anionic hydrophilic cream SR...	4061406	1
Apomix/Noweda	Anionic hydrophlic cream wit...	3B240214	1
Gehe/Apomix	Anionic hydrophilic cream SR...	4A030214	2
AMH Niemann / E	Anionic hydrophilic cream SR...	4A030214	1
phönix	Anionic hydrophilic cream SR...	4a130217	1
Apomix /Noweda	Anionic hydrophilic cream SR...	4A130217	1
Apomix/Noweda	Anionic hydrophilic cream SR...	4A270115	1
Apomix	Anionic hydrophlic cream wit...	4B020517	1
AMH Niemann / G	Anionic hydrophlic cream wit...	4B100314	1
APOMIX/Noweda	Anionic hydrophlic cream wit...	4B210316	1
PKH Halle	Anionic hydrophlic cream wit...	4B210316	1
apomix/Phönix	Anionic hydrophlic cream wit...	4B250315	1
PKH Halle / Phoenix	Anionic hydrophilic cream SR...	5A170214	1
apomix/ Noweda	Anionic hydrophilic cream SR...	5A200116	1
Allianc. Health. / Apomix	Anionic hydrophilic cream SR...	5A200217	3
Caelo	Anionic hydrophilic cream SR...	6061502	1
Apomix/Noweda	Anionic hydrophlic cream wit...	5B070414	1
All.Healthc. /Apomix	Anionic hydrophilic cream SR...	6A050215	1
All. Healthc. / Apomix	Anionic hydrophilic cream SR...	6050215	1
Anzag	Anionic hydrophilic cream SR...	6A260116	1
Apomix/Kehr	Anionic hydrophilic cream SR...	6A280217	1
Apomix/Noweda	Anionic hydrophlic cream wit...	6B040515	1
Apomix	Anionic hydrophlic cream wit...	6B040515	1
Apomix	Anionic hydrophlic cream wit...	6B04515	1
Noweda	Anionic hydrophlic cream wit...	6B040515	1
AMH Niemann / G	Anionic hydrophlic cream wit...	6B060514	1
apomix/Noweda	Anionic hydrophlic cream wit...	6B060514	1
AMH Niemann / G	Anionic hydrophlic cream wit...	6B190516	1
APOMIX/Noweda	Anionic hydrophilic cream SR...	7a090216	1
phönix	Anionic hydrophilic cream SR...	7a090216	1
Apomix	Anionic hydrophilic cream SR...	7A140317	1
Apomix/Phönix	Anionic hydrophilic cream SR...	7A200318	1
AMH Niemann / G	Anionic hydrophlic cream wit...	7B080615	1
apomix	Anionic hydrophlic cream wit...	7B210616	1
Apomix	Anionic hydrophilic cream SR...	8A090216	1
Apomix/Noweda	Anionic hydrophilic cream SR...	8A260314	1
Caelo	Anionic hydrophilic cream wi...	9771644	3
PKH Halle/Kehr	Anionic hydrophilic cream SR...	9A170216	1
AMH Niemann / G	Anionic hydrophilic cream SR...	9A210317	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix/Phönix	Anionic hydrophilic cream SR...	9A210317	1
AMW Niemann GmbH & Co. KG	Anionic hydrophlic cream wit...	9B150816	1
Apomix/Noweda	Anionic hydrophlic cream wit...	9B130715	1
Apomix	Anionic hydrophilic cream SR...	10A010414	1
Apomix	Anionic hydrophlic cream wit...	9B181017	1
Apomix/Noweda	Anionic hydrophlic cream wit...	9B220714	1
Caelo	Anionic hydrophilic cream DAB	108/AKD16	1
Apomix/Noweda	Anionic hydrophilic cream SR...	10A260318	1
AMH Niemann / G	Anionic hydrophilic cream SR...	10A050315	1
phönix	Anionic hydrophilic cream SR...	10a280317	1
Apomix/Phönix	Anionic hydrophilic cream SR...	10A280317	1
apomix/Holdermann	Anionic hydrophilic cream SR...	10A280317	1
AMH Niemann / G	Anionic hydrophlic cream wit...	10B060814	1
Apomix	Anionic hydrophlic cream wit...	10B060916	1
	Anionic hydrophilic cream SR...	110913	1
Pharmachem/Phönix	Anionic hydrophilic cream SR...	110913	1
Caelo	Anionic hydrophilic cream DAB	11170717	1
apomix/Phönix	Anionic hydrophilic cream SR...	11A010316	1
Apomix	Anionic hydrophilic cream SR...	11A010316	1
Caelo	Anionic hydrophilic cream SR...	11A010414	1
Raben-Apotheke Mildenau	Anionic hydrophilic cream SR...	11A010414	1
AMH Niemann / G	Anionic hydrophilic cream SR...	11A010414	1
Apomix/Sanacorp	Anionic hydrophilic cream SR...	11A190417	1
apomix/GEHE/28.05.14	Anionic hydrophilic cream SR...	11A010414	1
Caelo	Anionic hydrophilic cream DAB	1234	1
Caelo	Anionic hydrophilic cream DAB	12345	1
Caelo	Anionic hydrophilic cream DAB	123456	1
Caelo	Anionic hydrophilic cream DAB	123test	1
Aallianc. Healthc. / Ap...	Anionic hydrophilic cream SR...	12A030517	1
Apomix/Noweda	Anionic hydrophilic cream SR...	12A310315	1
Apomix	Anionic hydrophilic cream SR...	12A080414	1
AMH Niemann / G	Anionic hydrophilic cream SR...	12A310315	1
apomix/Phönix	Anionic hydrophilic cream SR...	12A3103215	1
apomix/Phönix	Anionic hydrophilic cream SR...	12A310315	1
Apomix/Noweda	Anionic hydrophlic cream wit...	12B212015	1
Apomix	Anionic hydrophlic cream wit...	12B220914	2
Apomix/Noweda	Anionic hydrophlic cream wit...	12B220914	1
Caelo	Anionic hydrophilic cream wi...	13067601	1
Pharmachem/Phönix	Anionic hydrophilic cream SR...	130814	1
PKH Halle	Anionic hydrophilic cream SR...	13A010415	1
Apomix/Noweda	Anionic hydrophilic cream SR...	13A080414	1
Apomix	Anionic hydrophilic cream SR...	13A030517	1
AMH Niemann / G	Anionic hydrophilic cream SR...	13A111016	1
AMH Niemann / K	Anionic hydrophilic cream SR...	13A220316	1
Apomix	Anionic hydrophlic cream wit...	13B141216	1
PKH Halle	Anionic hydrophlic cream wit...	13B141216	1
apomix/Phönix	Anionic hydrophlic cream wit...	13B161115	1
Caelo	Anionic hydrophilic cream SR...	14063303	2
Caelo	Anionic hydrophilic cream SR...	14041403	1
Caelo	Anionic hydrophilic cream SR...	14182101	2
Caelo	Anionic hydrophilic cream SR...	14244701	6
Caelo	Anionic hydrophilic cream SR...	14186501	2
AHD	Anionic hydrophilic cream SR...	14244702	1
AHD	Anionic hydrophilic cream SR...	14301501	2
Caelo	Anionic hydrophilic cream SR...	14301501	4
Caelo	Anionic hydrophilic cream wi...	14334501	19
Caelo	Anionic hydrophilic cream SR...	14336201	1
Caelo	Anionic hydrophilic cream SR...	14336203	3
Apomix	Anionic hydrophilic cream SR...	14A070514	3
phönix	Anionic hydrophilic cream SR...	14a150517	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
AMH Niemann / G	Anionic hydrophilic cream SR...	14A300316	3
Caelo	Anionic hydrophilic cream SR...	15014801	4
Caelo	Anionic hydrophilic cream SR...	15014804	1
Gehe	Anionic hydrophilic cream SR...	15014804	1
Caelo	Anionic hydrophilic cream SR...	15069403	1
Caelo	Anionic hydrophilic cream SR...	15069501	2
Caelo	Anionic hydrophilic cream SR...	15134602	2
Bombastus	Anionic hydrophilic cream SR...	150820UEA	1
Pharmachemie Pößneck	Anionic hydrophilic cream SR...	150716	1
Fagron	Anionic hydrophilic cream SR...	15172901	1
Caelo	Anionic hydrophilic cream SR...	15172901	5
Caelo	Anionic hydrophilic cream SR...	15205501	3
Caelo	Anionic hydrophilic cream SR...	15205502	2
Caelo	Anionic hydrophilic cream SR...	15268001	1
Caelo	Anionic hydrophilic cream DAB	15250201	1
Caelo	Anionic hydrophilic cream SR...	15268002	1
Caelo	Anionic hydrophilic cream DAB	15315601	3
Caelo	Anionic hydrophilic cream wi...	15272902	6
sanacorp	Anionic hydrophilic cream SR...	15319101	1
Caelo	Anionic hydrophilic cream SR...	15329501	7
Caelo	Anionic hydrophilic cream SR...	15329601	7
Noweda	Anionic hydrophilic cream SR...	15329601	1
Noweda	Anionic hydrophilic cream SR...	15386502	1
Caelo	Anionic hydrophilic cream SR...	15386502	3
Caelo	Anionic hydrophilic cream SR...	15386602	2
Caesar&Loretz GmbH	Anionic hydrophilic cream SR...	15386602	1
Caelo	Anionic hydrophilic cream DAB	15400701	7
Caelo	Anionic hydrophilic cream wi...	15405302	4
Caelo	Anionic hydrophilic cream DAB	15431701	9
Caelo	Anionic hydrophilic cream DAB	15428501	4
Apomix/Noweda	Anionic hydrophilic cream SR...	15A060416	1
Caesar & Loretz GmbH Hi...	Anionic hydrophilic cream DAB	15432001	2
Apomix/Noweda	Anionic hydrophilic cream SR...	15A120514	1
Apomix/Noweda	Anionic hydrophilic cream wit...	15B291214	1
PKH / Phoenix	Anionic hydrophilic cream SR...	15A210415	1
Caelo	Anionic hydrophilic cream SR...	160107502	1
Caelo	Anionic hydrophilic cream SR...	16019601	2
Caelo	Anionic hydrophilic cream SR...	16019801	5
Caesar&Loretz GmbH	Anionic hydrophilic cream SR...	16019601	1
Caelo	Anionic hydrophilic cream DAB	16046201	7
Caesar&Loretz GmbH,Phön...	Anionic hydrophilic cream SR...	16046501	1
Caelo	Anionic hydrophilic cream SR...	16046601	2
Caelo	Anionic hydrophilic cream DAB	16046901	10
Caelo	Anionic hydrophilic cream DAB	16046601	2
Euro OTC	Anionic hydrophilic cream DAB	16046701	1
Caelo	Anionic hydrophilic cream DAB	160470001	1
Caelo	Anionic hydrophilic cream DAB	16047001	10
Caelo	Anionic hydrophlic cream wit...	16047001	1
Caelo	Anionic hydrophilic cream DAB	16057201	5
Fishar	Anionic hydrophilic cream DAB	16057201	1
Caelo	Anionic hydrophilic cream DAB	16057301	9
Caelo	Anionic hydrophilic cream DAB	16057401	5
Caelo	Anionic hydrophilic cream DAB	16057501	4
Caelo	Anionic hydrophilic cream DAB	16057601	5
Caelo	Anionic hydrophilic cream DAB	16057701	6
Caelo	Anionic hydrophilic cream DAB	16057901	2
Caelo	Anionic hydrophilic cream DAB	16058001	5
Caelo	Anionic hydrophlic cream wit...	16058001	2
Caelo	Anionic hydrophilic cream wi...	16067301	5
Gehe	Anionic hydrophilic cream wi...	16067401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream wi...	16067401	7
Caelo	Anionic hydrophilic cream wi...	16067501	7
Caelo	Anionic hydrophilic cream wi...	16067601	23
Caelo/Phoenix	Anionic hydrophilic cream wi...	16067601	1
Caeanlo	Anionic hydrophilic cream wi...	16067601	1
Bombastus	Anionic hydrophilic cream wi...	16067701	2
Caelo	Anionic hydrophilic cream wi...	16067701	10
Caelo	Anionic hydrophilic cream wi...	16067801	6
Martinus-Apotheke	Anionic hydrophilic cream wi...	16067801	1
Caelo	Anionic hydrophilic cream wi...	16067901	3
Caelo	Anionic hydrophilic cream DAB	160914UEA	1
Caesar&lorenz,phönix,11...	Anionic hydrophilic cream SR...	16103101	3
Caelo	Anionic hydrophilic cream DAB	16103101	5
wepa/Sanaco	Anionic hydrophilic cream DAB	16103301	1
Caelo	Anionic hydrophilic cream DAB	16103301	8
Klenk	Anionic hydrophilic cream DAB	16103301	2
Caelo	Anionic hydrophilic cream DAB	16103401	8
Caelo	Anionic hydrophilic cream DAB	1610341	1
Caelo	Anionic hydrophilic cream DAB	16103501	3
Fagron	Anionic hydrophilic cream DAB	16103501	1
Cealo/AHD	Anionic hydrophilic cream SR...	16107501	1
Caelo	Anionic hydrophilic cream SR...	16107501	11
Caelo	Anionic hydrophilic cream SR...	16107502	5
Caesar&Loretz GmbH	Anionic hydrophilic cream SR...	16107502	1
Caelo	Anionic hydrophilic cream SR...	16107503	2
Caelo	Anionic hydrophilic cream SR...	16114102	1
Caelo	Anionic hydrophilic cream DAB	161101UEA/16156401	1
Caelo	Anionic hydrophilic cream DAB	16150701	4
Caelo	Anionic hydrophilic cream SR...	16126701	1
Ichthyol/Phönix	Anionic hydrophilic cream DAB	16150702	1
Caelo	Anionic hydrophilic cream DAB	16150702	8
Caelo	Anionic hydrophilic cream DAB	16150703	1
Caelo	Anionic hydrophilic cream DAB	16150801	11
Euro OTC	Anionic hydrophilic cream DAB	16150801	1
Caelo	Anionic hydrophilic cream SR...	16150901	3
Caelo	Anionic hydrophilic cream DAB	16150901	5
Caelo	Anionic hydrophilic cream DAB	16152501	3
Caelo	Anionic hydrophilic cream DAB	16152701	10
Sanacorp	Anionic hydrophilic cream DAB	16152701	1
Caelo	Anionic hydrophilic cream DAB	16152801	7
Caelo	Anionic hydrophilic cream DAB	16152901	1
Caelo	Anionic hydrophilic cream DAB	16156001	6
Caelo	Anionic hydrophilic cream DAB	16156101	8
Noweda	Anionic hydrophilic cream DAB	16156401	2
Caelo	Anionic hydrophilic cream DAB	16156401	10
Caelo	Anionic hydrophilic cream DAB	16158501	1
Alte Apotheke	Anionic hydrophilic cream DAB	16158501	1
Caelo	Anionic hydrophilic cream DAB	16158701	13
phönix	Anionic hydrophilic cream SR...	16158802	1
Caelo	Anionic hydrophilic cream DAB	16158803	4
Caelo	Anionic hydrophilic cream DAB	16158901	11
Caesar & Loretz/Noweda	Anionic hydrophilic cream DAB	16158901	1
Caelo	Anionic hydrophilic cream DAB	16179302	1
Caelo	Anionic hydrophilic cream DAB	16179601	8
Gehe	Anionic hydrophilic cream DAB	16179601	1
Caelo	Anionic hydrophilic cream SR...	16179601	1
Caelo	Anionic hydrophilic cream DAB	16179702	5
Caelo	Anionic hydrophilic cream DAB	16179801	9
Caelo	Anionic hydrophilic cream DAB	16179901	2
Caelo	Anionic hydrophilic cream DAB	16180001	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp 28.11 Ek.5,85	Anionic hydrophilic cream DAB	16180001	1
Caelo	Anionic hydrophilic cream DAB	16180101	2
Phoenix	Anionic hydrophilic cream SR...	16188903	1
Caelo	Anionic hydrophilic cream wi...	16191501	2
Gehe	Anionic hydrophilic cream wi...	16191501	1
Caelo	Anionic hydrophilic cream wi...	16191601	17
Noweda	Anionic hydrophilic cream wi...	16191601	1
Caelo	Anionic hydrophilic cream wi...	16191701	5
Caelo	Anionic hydrophilic cream wi...	16191801	4
Caelo	Anionic hydrophilic cream SR...	162115501	1
Caelo	Anionic hydrophilic cream SR...	16215501	6
Caelo	Anionic hydrophilic cream wi...	16191901	4
Caelo	Anionic hydrophlic cream wit...	16215501	1
Caeanto	Anionic hydrophilic cream SR...	16215501	2
Caelo	Anionic hydrophilic cream SR...	16215601	3
Caelo	Anionic hydrophilic cream SR...	16215701	1
Caelo	Anionic hydrophilic cream DAB	16220201	1
Caelo	Anionic hydrophilic cream DAB	16249901	1
Caelo	Anionic hydrophilic cream DAB	16249902	1
Celo	Anionic hydrophilic cream DAB	16250001	4
Bombastus	Anionic hydrophilic cream DAB	16250001	1
Caelo	Anionic hydrophilic cream DAB	16250001	13
Hochhaus-Apotheke	Anionic hydrophilic cream DAB	16250001	2
gepepharm/Gehe	Anionic hydrophilic cream DAB	16250001	1
Caelo	Anionic hydrophilic cream DAB	16250201	10
Caelo	Anionic hydrophilic cream SR...	16258601	7
Caelo	Anionic hydrophilic cream DAB	16258601	2
Caelo	Anionic hydrophilic cream SR...	16258701	6
Caelo	Anionic hydrophilic cream DAB	16259001	5
Noweda	Anionic hydrophilic cream DAB	16259001	1
Caelo	Anionic hydrophilic cream DAB	1625901	1
Caelo	Anionic hydrophilic cream DAB	16259301	13
Calo/Sanacorp	Anionic hydrophilic cream DAB	16259302	1
Caelo	Anionic hydrophilic cream DAB	16259302	1
Caelo	Anionic hydrophilic cream DAB	16259401	8
Caelo	Anionic hydrophilic cream DAB	16259601	5
Caelo	Anionic hydrophilic cream wi...	1625961	1
Caelo	Anionic hydrophilic cream DAB	16259701	10
20.12.2016 Sanacorp 6,00EUR	Anionic hydrophilic cream DAB	16259701	1
Caelo	Anionic hydrophilic cream DAB	16259801	5
Caelo	Anionic hydrophilic cream SR...	16259801	1
Caelo	Anionic hydrophilic cream DAB	16259901	9
Caelo	Anionic hydrophilic cream DAB	16260001	5
Caelo	Anionic hydrophilic cream wi...	16260001	1
Caelo	Anionic hydrophilic cream DAB	16264301	8
Caelo	Anionic hydrophilic cream DAB	16264401	1
Caelo	Anionic hydrophilic cream DAB	16264501	8
Sanacorp	Anionic hydrophilic cream DAB	16264501	1
Caelo	Anionic hydrophilic cream DAB	16264601	3
Caelo	Anionic hydrophilic cream DAB	16264701	5
Caelo	Anionic hydrophilic cream SR...	16289402	1
Caelo	Anionic hydrophilic cream wi...	16290101	8
Caelo	Anionic hydrophilic cream wi...	16290201	1
Caelo	Anionic hydrophilic cream wi...	16290301	14
Caelo	Anionic hydrophilic cream wi...	16290302	4
Caelo	Anionic hydrophilic cream wi...	16290401	5
Caelo	Anionic hydrophilic cream wi...	16290501	7
Caelo	Anionic hydrophilic cream DAB	163170001	1
Caelo	Anionic hydrophilic cream DAB	16317001	3
Sanacorp WE 23.02.17 -1...	Anionic hydrophilic cream wi...	16290501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	16317201	3
Noweda	Anionic hydrophilic cream DAB	16317201	1
Caelo	Anionic hydrophilic cream DAB	16322401	6
Finteler Apotheke	Anionic hydrophilic cream DAB	16322402	1
Caelo	Anionic hydrophilic cream DAB	16322501	2
Caelo	Anionic hydrophilic cream DAB	16322801	10
Caesar&LoretzGmbH,11,03...	Anionic hydrophilic cream SR...	16322701	1
Phoenix	Anionic hydrophilic cream DAB	16322801	1
Caelo	Anionic hydrophilic cream DAB	16328001	5
AHD 09.03.2017 EK: 5,94	Anionic hydrophilic cream DAB	16328001	1
Allianz Healthcare	Anionic hydrophilic cream DAB	16328001	1
Caelo	Anionic hydrophilic cream DAB	16328201	6
Caelo	Anionic hydrophilic cream DAB	16328101	1
Caelo	Anionic hydrophilic cream DAB	16328401	7
Allianz Healthcare	Anionic hydrophilic cream DAB	16328401	1
Caelo	Anionic hydrophilic cream SR...	16332001	5
Caelo	Anionic hydrophilic cream SR...	16332101	2
Caelo	Anionic hydrophilic cream DAB	1659801	1
AMH Niemann / E	Anionic hydrophilic cream SR...	16A050515	1
Gehe/Apomix	Anionic hydrophilic cream SR...	16A130514	1
apomix/ Noweda	Anionic hydrophilic cream SR...	16A060416	1
PKH/Gehe	Anionic hydrophilic cream SR...	16A060416	1
Bombastus	Anionic hydrophilic cream SR...	17000162	3
Gehe	Anionic hydrophilic cream SR...	17000163	1
Bombastus	Anionic hydrophilic cream SR...	17000163	1
Bombastus	Anionic hydrophilic cream SR...	17000166	3
Bombastus	Anionic hydrophilic cream SR...	17000164	2
Bombastus	Anionic hydrophilic cream DAB	17001302	3
Bombastus	Anionic hydrophilic cream DAB	17001305	1
Bombastus	Anionic hydrophilic cream DAB	17001313	1
Bombastus	Anionic hydrophilic cream SR...	17001699	3
Bombastus	Anionic hydrophilic cream DAB	17001307	5
Caelo	Anionic hydrophilic cream SR...	17008901	2
Bombastus	Anionic hydrophilic cream SR...	17006166	1
Bombastus	Anionic hydrophilic cream SR...	17001702	1
Bombastus	Anionic hydrophilic cream SR...	170103CAH/303735	1
Caelo	Anionic hydrophilic cream wi...	17029901	7
Phoenix	Anionic hydrophilic cream wi...	17029901	1
Caelo	Anionic hydrophilic cream wi...	17029902	3
Caelo	Anionic hydrophilic cream wi...	17030001	4
Caelo	Anionic hydrophilic cream wi...	17030101	2
Sanacorp	Anionic hydrophilic cream wi...	17030202	2
Caelo	Anionic hydrophilic cream wi...	17030202	2
Caelo	Anionic hydrophilic cream wi...	17030301	15
Sanacorp	Anionic hydrophilic cream wi...	17030301	1
Allianz Healthcare	Anionic hydrophilic cream wi...	17030301	1
Bombastus	Anionic hydrophilic cream SR...	170313CAH/303736	1
Caelo	Anionic hydrophilic cream DAB	170411AHC/16264501	1
Caelo	Anionic hydrophilic cream SR...	17041701	1
Caelo	Anionic hydrophilic cream SR...	17041901	2
Caelo	Anionic hydrophilic cream SR...	17041902	3
Bombastus	Anionic hydrophilic cream DAB	17045101	1
Caelo	Anionic hydrophilic cream DAB	17045101	2
Caelo	Anionic hydrophilic cream DAB	17045201	2
Caelo	Anionic hydrophilic cream DAB	17045301	7
Kurpark Apo Defektur	Anionic hydrophilic cream DAB	17045301	1
Caesar&Loretz GmbH,26.0...	Anionic hydrophilic cream SR...	17045401	1
Caelo	Anionic hydrophilic cream DAB	17045401	4
Caelo	Anionic hydrophilic cream DAB	17045501	6
Fagron	Anionic hydrophilic cream DAB	17045501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	170503AHC/16264501	1
Caelo	Anionic hydrophilic cream DAB	170612UEA/16264501	1
Caelo	Anionic hydrophilic cream DAB	170817UEA/17088801	1
Caelo	Anionic hydrophilic cream DAB	17083801	1
Caelo	Anionic hydrophilic cream DAB	17083802	3
Caelo	Anionic hydrophilic cream DAB	17084101	12
Caelo	Anionic hydrophilic cream DAB	17084001	1
Noweda	Anionic hydrophilic cream DAB	17084101	2
Sanacorp	Anionic hydrophilic cream DAB	17084201	1
Caelo	Anionic hydrophilic cream DAB	17088701	1
apomix/Noweda	Anionic hydrophilic cream DAB	170841901	1
Caelo	Anionic hydrophilic cream SR...	17084201	1
Caelo	Anionic hydrophilic cream DAB	17084201	3
Caelo	Anionic hydrophilic cream DAB	17088801	1
Caesar&Loretz GmbH,PHön...	Anionic hydrophilic cream SR...	17088801	1
Caelo	Anionic hydrophilic cream wi...	17088901	1
Caesar & Loretz/Noweda	Anionic hydrophilic cream DAB	17088901	1
Caelo	Anionic hydrophilic cream DAB	17088901	2
Finteler Apotheke	Anionic hydrophilic cream DAB	17088901	1
Caelo	Anionic hydrophilic cream DAB	170919UEA/17143001	1
Caelo	Anionic hydrophilic cream SR...	17101501	2
Caelo	Anionic hydrophilic cream SR...	17101502	4
Caelo	Anionic hydrophilic cream SR...	17101504	1
Caelo	Anionic hydrophilic cream wi...	17102001	3
Caelo	Anionic hydrophilic cream wi...	17102101	2
Caelo	Anionic hydrophilic cream wi...	17102201	4
Caelo	Anionic hydrophilic cream wi...	17102301	3
Caelo	Anionic hydrophilic cream wi...	17102301	2
Caelo	Anionic hydrophilic cream wi...	17102401	3
Caelo	Anionic hydrophilic cream DAB	171107UEA/17143001	1
Caelo	Anionic hydrophilic cream DAB	171229UEA/17123001	1
Caelo	Anionic hydrophilic cream DAB	17125802	1
Caelo	Anionic hydrophilic cream SR...	17125802	1
Caelo	Anionic hydrophilic cream DAB	17141901	5
Caelo	Anionic hydrophilic cream DAB	171420101	1
Caelo	Anionic hydrophilic cream DAB	17142501	1
Caelo	Anionic hydrophilic cream DAB	17142601	2
Caesar & Loretz GmbH, P...	Anionic hydrophilic cream SR...	17142801	1
Caelo	Anionic hydrophilic cream DAB	17142801	1
Caelo	Anionic hydrophilic cream DAB	17142901	2
Caelo	Anionic hydrophilic cream DAB	17143001	2
Caelo	Anionic hydrophilic cream DAB	17143101	1
Noweda	Anionic hydrophilic cream DAB	17143101	1
Noweda	Anionic hydrophilic cream wit...	17143001	3
Sanacorp	Anionic hydrophilic cream DAB	17143101	1
Caelo	Anionic hydrophilic cream SR...	17143201	1
Sanacorp	Anionic hydrophilic cream DAB	17143201	1
Caelo	Anionic hydrophilic cream DAB	17143301	1
Caelo	Anionic hydrophilic cream DAB	17143601	1
Caelo	Anionic hydrophilic cream DAB	17143701	3
Caelo	Anionic hydrophilic cream DAB	17143401	1
Caelo	Anionic hydrophilic cream DAB	17143801	1
Caelo	Anionic hydrophilic cream DAB	17143802	4
Caelo	Anionic hydrophilic cream wi...	17156501	3
Caelo	Anionic hydrophilic cream SR...	17153802	1
Caelo	Anionic hydrophilic cream wi...	17156601	1
Caelo	Anionic hydrophilic cream SR...	17169301	2
Noweda	Anionic hydrophilic cream wi...	17156601	1
Caelo	Anionic hydrophilic cream wi...	17185401	1
Caelo	Anionic hydrophilic cream SR...	17184201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	17200701	1
Caelo	Anionic hydrophilic cream DAB	17200703	1
Caelo	Anionic hydrophilic cream DAB	17201102	1
Caelo	Anionic hydrophilic cream DAB	17201201	2
Caelo	Anionic hydrophilic cream DAB	17201301	5
Caelo	Anionic hydrophilic cream DAB	17201401	5
Caelo	Anionic hydrophilic cream DAB	17201601	1
Caelo	Anionic hydrophilic cream DAB	17201701	1
Caelo	Anionic hydrophilic cream DAB	17201501	3
Caelo	Anionic hydrophilic cream DAB	17201801	1
Caesar & Loretz/AHD	Anionic hydrophilic cream DAB	17202001	1
Caelo	Anionic hydrophilic cream DAB	17202001	2
Caelo	Anionic hydrophilic cream DAB	17202101	1
Caelo	Anionic hydrophilic cream DAB	17202301	2
Caelo	Anionic hydrophilic cream DAB	17208001	1
Caelo	Anionic hydrophilic cream SR...	17211701	1
Caelo	Anionic hydrophilic cream SR...	17211702	2
Caelo	Anionic hydrophilic cream DAB	1724601	1
Caelo	Anionic hydrophilic cream SR...	17256901	3
Caelo	Anionic hydrophilic cream DAB	17248701	1
Caelo	Anionic hydrophilic cream wi...	17261901	1
Caelo	Anionic hydrophilic cream wi...	17262002	1
Caelo	Anionic hydrophilic cream wi...	1726202	1
Caelo	Anionic hydrophilic cream wi...	17262101	6
Fagron	Anionic hydrophilic cream wi...	17262101	1
Caelo	Anionic hydrophilic cream wi...	17262201	1
Caelo	Anionic hydrophilic cream DAB	17270001	5
Caelo	Anionic hydrophilic cream DAB	17270301	1
Caelo	Anionic hydrophilic cream DAB	17270401	2
Caelo	Anionic hydrophilic cream wi...	17270501	1
Caelo	Anionic hydrophilic cream DAB	17270501	3
Caelo	Anionic hydrophilic cream SR...	17294401	1
Lamotte	Anionic hydrophilic cream DAB	17299101	1
Caelo	Anionic hydrophilic cream SR...	17294501	1
Caelo	Anionic hydrophilic cream DAB	17299101	3
Fagron	Anionic hydrophilic cream DAB	17299201	1
Caelo	Anionic hydrophilic cream DAB	17299201	1
Caesar & Loretz/ Noweda	Anionic hydrophilic cream DAB	17299201	1
Caelo	Anionic hydrophilic cream DAB	17299801	1
Caelo	Anionic hydrophilic cream DAB	17299501	1
Caelo	Anionic hydrophilic cream SR...	17299701	1
Caelo	Anionic hydrophilic cream wi...	17335401	3
Caelo	Anionic hydrophilic cream wi...	17335501	1
Caelo	Anionic hydrophilic cream wi...	17335701	1
Caelo	Anionic hydrophilic cream wi...	17335601	1
Apomix/AHD	Anionic hydrophilic cream SR...	17A030516	1
Apomix/Sanacorp	Anionic hydrophilic cream SR...	17A030516	1
PKH Halle / Phoenix	Anionic hydrophilic cream SR...	17A200514	1
Noweda	Anionic hydrophilic cream SR...	17A120515	1
Bombastus	Anionic hydrophilic cream DAB	18000257	3
Bombastus	Anionic hydrophilic cream SR...	18001475	3
Bombastus	Anionic hydrophilic cream SR...	18001482	1
Bombastus	Anionic hydrophilic cream SR...	18001474	1
Bombastus	Anionic hydrophilic cream DAB	18001549	1
EuRho/Noweda	Anionic hydrophilic cream SR...	18001681	1
Bombastus	Anionic hydrophilic cream DAB	18002356	1
Caelo	Anionic hydrophilic cream SR...	18041104	5
Caelo	Anionic hydrophilic cream DAB	180220UEA/17200801	1
Caelo	Anionic hydrophilic cream SR...	18041105	1
Caelo	Anionic hydrophilic cream wi...	18041501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream SR...	18067801	2
Allianc. Healthc. /Apomix	Anionic hydrophilic cream SR...	18A040516	2
Apomix/Noweda	Anionic hydrophilic cream SR...	18A040516	1
Apomix	Anionic hydrophilic cream SR...	18A040516	1
Phönix	Anionic hydrophilic cream SR...	18A140617	1
Apomix	Anionic hydrophilic cream SR...	18A140617	1
apomix/Noweda	Anionic hydrophilic cream SR...	18A260514	2
Apomix/Noweda	Anionic hydrophilic cream SR...	18A260514	1
All..Healthc./Apomix	Anionic hydrophilic cream SR...	18A260514	1
Alliance Healthcare/Apo...	Anionic hydrophilic cream SR...	19A030614	1
Raben-Apotheke Mildenau	Anionic hydrophilic cream SR...	19A030614	1
Gehe/apomix/23.09.14	Anionic hydrophilic cream SR...	19A030614	1
Apomix/Gehe	Anionic hydrophilic cream SR...	19A260515	2
apomix/Gehe	Anionic hydrophilic cream SR...	19A030614	1
Apomix/Kehr	Anionic hydrophilic cream SR...	19A180516	1
Apomix/Noweda	Anionic hydrophilic cream SR...	20A020615	2
AMH Niemann / E	Anionic hydrophilic cream SR...	20A020615	1
Apomix	Anionic hydrophilic cream SR...	21A080616	1
Noweda/apomix	Anionic hydrophilic cream SR...	21A090615	1
Apomix/Phönix	Anionic hydrophilic cream SR...	21A180614	1
Apomix/Noweda	Anionic hydrophilic cream SR...	21A180614	1
Apomix	Anionic hydrophilic cream SR...	21A220617	1
Bombastus	Anionic hydrophilic cream SR...	2215M-03327	1
Apomix/Phönix	Anionic hydrophilic cream SR...	22A020817	1
Apomix/Phönix	Anionic hydrophilic cream SR...	22A040717	1
Apomix	Anionic hydrophilic cream SR...	22A040717	1
Noweda	Anionic hydrophilic cream SR...	22A040717	1
AMH Niemann / K	Anionic hydrophilic cream SR...	22A090616	1
Apomix	Anionic hydrophilic cream SR...	22A090616	1
Allianc. Health. /Apomix	Anionic hydrophilic cream SR...	22A090616	2
Caelo	Anionic hydrophilic cream SR...	22A160615	2
Apomix	Anionic hydrophilic cream SR...	23A130616	2
Caelo	Anionic hydrophilic cream SR...	23121411	1
Allianc. Health. /Apomix	Anionic hydrophilic cream SR...	23A130616	2
apomix	Anionic hydrophilic cream SR...	23A130616	1
Gehe/02.09.14/apomix	Anionic hydrophilic cream SR...	23A240614	1
Pharmachemie	Anionic hydrophilic cream SR...	240915	1
apomix/AHD	Anionic hydrophilic cream SR...	24A280616	1
Noweda/Apomix	Anionic hydrophilic cream SR...	25A290616	1
Neidhardt	Anionic hydrophilic cream SR...	25A140714	1
PKH Halle/Kehr	Anionic hydrophilic cream SR...	24A280616	1
Euro OTC	Anionic hydrophilic cream DAB	260416v2	1
apomix/Phönix	Anionic hydrophilic cream SR...	25A290616	1
Pharmachem/Phoenix	Anionic hydrophilic cream SR...	260712	3
Gehe/apomix/17.12.2014	Anionic hydrophilic cream SR...	26A060814	1
Apomix/Phönix	Anionic hydrophilic cream SR...	26A130716	2
apomix/Holdermann	Anionic hydrophilic cream SR...	26A130716	1
Bombastus	Anionic hydrophilic cream SR...	280568	2
Caelo	Anionic hydrophilic cream SR...	28072014A	1
Bombastus	Anionic hydrophilic cream SR...	281738	7
Bombastus	Anionic hydrophilic cream SR...	281740	5
Bombastus	Anionic hydrophilic cream SR...	281767	6
Bombastu-Werke AG/Gehe	Anionic hydrophilic cream SR...	281866	2
Bombastus	Anionic hydrophilic cream SR...	281866	9
Bombastus	Anionic hydrophilic cream SR...	283279	2
Bombastus	Anionic hydrophilic cream SR...	283283	3
Bombastus	Anionic hydrophilic cream SR...	283286	2
Bombastus	Anionic hydrophilic cream SR...	283285	4
Bombastus	Anionic hydrophilic cream SR...	284078	1
Bombastu-Werke AG/Gehe	Anionic hydrophilic cream SR...	283286	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream SR...	284079	1
Bombastus	Anionic hydrophilic cream SR...	284080	3
Bombastus	Anionic hydrophilic cream SR...	284081	5
Bombastus	Anionic hydrophilic cream SR...	284082	4
Bombastus	Anionic hydrophilic cream SR...	284083	4
Bombastus	Anionic hydrophilic cream SR...	285890	2
Phönix	Anionic hydrophilic cream SR...	284083	1
Bombastus	Anionic hydrophilic cream SR...	285829	1
Bombastus	Anionic hydrophilic cream SR...	285893	2
Bombastus	Anionic hydrophilic cream SR...	285892	1
Bombastus	Anionic hydrophilic cream SR...	285895	3
Bombastus	Anionic hydrophilic cream SR...	285894	2
Bombastus	Anionic hydrophilic cream SR...	285896	3
Bombastus	Anionic hydrophilic cream SR...	285897	3
Bombastus	Anionic hydrophilic cream SR...	285898	10
Bombastus	Anionic hydrophilic cream SR...	285899	5
anzag	Anionic hydrophilic cream SR...	285899	1
Cealo	Anionic hydrophilic cream SR...	285899	1
Euro OTC	Anionic hydrophilic cream SR...	286838	2
Bombastus	Anionic hydrophilic cream SR...	286972	1
Bombastus	Anionic hydrophilic cream SR...	286971	1
Bombastus	Anionic hydrophilic cream SR...	286973	1
Gehe	Anionic hydrophilic cream SR...	286974	1
Euro OTC	Anionic hydrophilic cream SR...	286975	1
Bombastus	Anionic hydrophilic cream SR...	286975	3
Caelo	Anionic hydrophilic cream SR...	286976	2
Bombastus	Anionic hydrophilic cream SR...	286977	2
Anzag	Anionic hydrophilic cream SR...	286979	1
Bombastus	Anionic hydrophilic cream SR...	286979	13
Bombastu-Werke AG/Gehe	Anionic hydrophilic cream SR...	286979	1
Caelo	Anionic hydrophilic cream SR...	286979	1
Bombastus	Anionic hydrophilic cream SR...	286980	9
Bombastus	Anionic hydrophilic cream SR...	288227	1
Bombastus	Anionic hydrophilic cream SR...	287717	1
Bombastus	Anionic hydrophilic cream SR...	288230	2
Bombastus	Anionic hydrophilic cream SR...	288231	4
Bombastus	Anionic hydrophilic cream SR...	288233	3
Bombastus	Anionic hydrophilic cream SR...	288234	2
Bombastus	Anionic hydrophilic cream SR...	288235	3
Caelo	Anionic hydrophilic cream SR...	288235	1
Bombastus	Anionic hydrophilic cream SR...	288307	6
Caelo	Anionic hydrophilic cream SR...	288307	2
Bombastus	Anionic hydrophilic cream SR...	288308	6
Bomastus	Anionic hydrophilic cream SR...	288308	1
Caelo	Anionic hydrophilic cream SR...	288310	1
Bombastus	Anionic hydrophilic cream SR...	288310	5
Bombastus	Anionic hydrophilic cream SR...	288882	1
Bombastus	Anionic hydrophilic cream SR...	289325	2
Bombastus	Anionic hydrophilic cream SR...	289326	2
Bombastus	Anionic hydrophilic cream SR...	289323	1
Bombastus	Anionic hydrophilic cream SR...	289322	1
Bombastus	Anionic hydrophilic cream SR...	289327	4
Bombastus	Anionic hydrophilic cream SR...	289328	2
Bombastus	Anionic hydrophilic cream SR...	289330	2
Bombastus	Anionic hydrophilic cream SR...	289331	4
Caelo	Anionic hydrophilic cream SR...	289331	1
Bombastus	Anionic hydrophilic cream SR...	289332	6
Caelo	Anionic hydrophilic cream SR...	289332	1
apomix/Phönix	Anionic hydrophilic cream SR...	28A040815	2
Apomix/Gehe	Anionic hydrophilic cream SR...	28A040815	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Anionic hydrophilic cream SR...	28A190814	1
apomix	Anionic hydrophilic cream SR...	28A190814	1
Bombastus	Anionic hydrophilic cream SR...	290522	1
Apomix	Anionic hydrophilic cream SR...	28A260716	2
Bombastus	Anionic hydrophilic cream SR...	290523	2
Bombastus	Anionic hydrophilic cream SR...	290525	1
Bombastus	Anionic hydrophilic cream SR...	290528	2
Bombastus	Anionic hydrophilic cream SR...	290524	1
Bombastus	Anionic hydrophilic cream SR...	290527	3
Dr. Wolff	Anionic hydrophilic cream SR...	290528	1
Bombastus	Anionic hydrophilic cream SR...	290530	7
Sanacorp	Anionic hydrophilic cream SR...	290530	1
Bombastus	Anionic hydrophilic cream SR...	291522	3
Bombastus	Anionic hydrophilic cream SR...	291523	7
Bombastus	Anionic hydrophilic cream SR...	291524	3
Bombastus	Anionic hydrophilic cream SR...	291525	3
Bombastus	Anionic hydrophilic cream SR...	291526	2
Bombastus	Anionic hydrophilic cream SR...	291527	3
Bombastus	Anionic hydrophilic cream SR...	291528	8
Bombastus	Anionic hydrophilic cream SR...	291529	7
Fagron	Anionic hydrophilic cream SR...	291528	1
Fagron	Anionic hydrophilic cream SR...	291529	1
Bombastus	Anionic hydrophilic cream SR...	292975	2
Bombastus	Anionic hydrophilic cream SR...	292976	5
Bombastus	Anionic hydrophilic cream SR...	292977	5
Bombastus	Anionic hydrophilic cream SR...	292978	4
Bombastus	Anionic hydrophilic cream SR...	292979	3
Bombastus	Anionic hydrophilic cream SR...	293102	2
Caelo	Anionic hydrophilic cream SR...	293103	1
Bombastus	Anionic hydrophilic cream SR...	293103	2
Caelo	Anionic hydrophilic cream SR...	293104	2
Bombastus	Anionic hydrophilic cream SR...	293104	5
Bombast/Phoenix	Anionic hydrophilic cream SR...	293104	1
Bombastus	Anionic hydrophilic cream SR...	293133	1
Bombastus	Anionic hydrophilic cream SR...	293134	2
Bombastus	Anionic hydrophilic cream SR...	293135	8
Caelo	Anionic hydrophilic cream SR...	293135	1
Bombastus	Anionic hydrophilic cream SR...	295165	1
Bombastus	Anionic hydrophilic cream SR...	295166	6
Bombastus	Anionic hydrophilic cream SR...	295167	5
Bombastus	Anionic hydrophilic cream SR...	295169	4
Bombastus	Anionic hydrophilic cream SR...	295170	2
Bombastus	Anionic hydrophilic cream SR...	295171	3
Bombastus	Anionic hydrophilic cream SR...	295232	8
Anzag	Anionic hydrophilic cream SR...	295232	1
Caelo	Anionic hydrophilic cream SR...	295232	1
Bombastus	Anionic hydrophilic cream SR...	295233	4
Bombastus	Anionic hydrophilic cream SR...	296406	1
Caelo	Anionic hydrophilic cream SR...	296408	1
Bombastus	Anionic hydrophilic cream SR...	296409	1
Euro OTC	Anionic hydrophilic cream SR...	296410	2
Bombastus	Anionic hydrophilic cream SR...	296411	2
Bombastus	Anionic hydrophilic cream SR...		1
Bombastus	Anionic hydrophilic cream SR...	296412	3
Bombastus	Anionic hydrophilic cream SR...	296413	2
Bombastus	Anionic hydrophilic cream SR...	296414	9
Bombastus	Anionic hydrophilic cream SR...	296415	3
Bombastus	Anionic hydrophilic cream SR...	297309	1
Bombastus	Anionic hydrophilic cream DAB	297667	1
Bombastus	Anionic hydrophilic cream SR...	297817	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream SR...	297819	1
Bombastus	Anionic hydrophilic cream SR...	297820	2
Bombastus	Anionic hydrophilic cream SR...	297821	5
Bombastus	Anionic hydrophilic cream SR...	297822	4
Bombastus	Anionic hydrophilic cream SR...	297823	5
Bombastus	Anionic hydrophilic cream SR...	297908	2
Bombastus	Anionic hydrophilic cream SR...	297910	3
Bombastus	Anionic hydrophilic cream SR...	297909	1
Bombastus	Anionic hydrophilic cream SR...	297911	3
Anzag	Anionic hydrophilic cream DAB	297911	2
Caelo	Anionic hydrophilic cream SR...	297912	1
Bombastus	Anionic hydrophilic cream SR...	297912	3
Bombastus	Anionic hydrophilic cream SR...	299253	4
Bombastus	Anionic hydrophilic cream SR...	297913	8
Bombastus	Anionic hydrophilic cream SR...	299254	3
Bombastus	Anionic hydrophilic cream SR...	299255	4
Bombastus	Anionic hydrophilic cream SR...	299256	2
Bombastus	Anionic hydrophilic cream DAB	299256	1
Apomix	Anionic hydrophilic cream SR...	29A020914	1
AMH Niemann / G	Anionic hydrophilic cream SR...	29A020914	1
apomix/Noweda	Anionic hydrophilic cream SR...	29A070917	1
Apomix /Noweda	Anionic hydrophilic cream SR...	29a160816	1
Apomix/Noweda	Anionic hydrophilic cream SR...	29A070917	1
Bombastus	Anionic hydrophilic cream SR...	2x297819	1
Caelo	Anionic hydrophilic cream SR...	30011503	1
Bombastus	Anionic hydrophilic cream DAB	300348	3
Bombastus	Anionic hydrophilic cream DAB	300349	2
Bombastus	Anionic hydrophilic cream DAB	300350	3
Bombastus	Anionic hydrophilic cream DAB	300353	3
Bombastus	Anionic hydrophilic cream DAB	300352	1
Bombastus	Anionic hydrophilic cream DAB	300351	4
Bombastus	Anionic hydrophilic cream DAB	300355	2
Bombastus	Anionic hydrophilic cream DAB	300357	2
Bombastus	Anionic hydrophilic cream DAB	300356	3
Bombastus	Anionic hydrophilic cream DAB	300358	5
Fagron	Anionic hydrophilic cream DAB	300358	1
Bombastus	Anionic hydrophilic cream DAB	300361	1
Bombastus	Anionic hydrophilic cream DAB	300363	1
Bombastus	Anionic hydrophilic cream DAB	300365	6
Caelo	Anionic hydrophilic cream DAB	300365	1
Bombastus	Anionic hydrophilic cream SR...	300560	1
Bombastus	Anionic hydrophilic cream SR...	300561	3
Bombastus	Anionic hydrophilic cream SR...	300562	2
Bombastus	Anionic hydrophilic cream SR...	300563	3
Euro OTC	Anionic hydrophilic cream SR...	300563	1
Bombastus	Anionic hydrophilic cream SR...	300564	15
Bombastus	Anionic hydrophilic cream SR...	300946	1
Bombastus	Anionic hydrophilic cream SR...	300947	1
Bombastus	Anionic hydrophilic cream SR...	300950	4
Bombastus	Anionic hydrophilic cream SR...	300951	4
Bombastus	Anionic hydrophilic cream SR...	300952	1
Bombastus	Anionic hydrophilic cream SR...	300953	3
Bombastus	Anionic hydrophilic cream SR...	300955	2
Caelo	Anionic hydrophilic cream SR...	300956	1
Bombastus	Anionic hydrophilic cream SR...	300956	9
Bombastus	Anionic hydrophilic cream SR...	300957	8
Caelo	Anionic hydrophilic cream SR...	300957	1
Bombastus	Anionic hydrophilic cream DAB	301267	4
Bombastus	Anionic hydrophilic cream DAB	301268	1
Bombastus	Anionic hydrophilic cream DAB	301269	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	301270	2
Bombastus	Anionic hydrophilic cream DAB	301271	3
Bombastus	Anionic hydrophilic cream DAB	301272	2
Bombastus	Anionic hydrophilic cream DAB	301349	3
Bombastus	Anionic hydrophilic cream DAB	301350	6
Bombastus	Anionic hydrophilic cream DAB	301351	3
Bombastus	Anionic hydrophilic cream DAB	301353	4
Bombastus	Anionic hydrophilic cream SR...	301352	1
Bombastus	Anionic hydrophilic cream DAB	301356	3
Bombastus	Anionic hydrophilic cream DAB	301357	2
Bombastus	Anionic hydrophilic cream DAB	301358	1
hochhaus-Apotheke	Anionic hydrophilic cream SR...	3016E-03195	1
Caelo	Anionic hydrophilic cream SR...	302196	2
Bombastus	Anionic hydrophilic cream SR...	302190	1
Bombastus	Anionic hydrophilic cream SR...	302194	4
Bombastus	Anionic hydrophilic cream SR...	302192	1
Bombastus	Anionic hydrophilic cream SR...	302191	1
Bombastus	Anionic hydrophilic cream SR...	302195	3
Bombastus	Anionic hydrophilic cream SR...	302196	2
Euro OTC	Anionic hydrophilic cream SR...	302199	1
Bombastus	Anionic hydrophilic cream SR...	302199	6
Caelo	Anionic hydrophilic cream SR...	302199	1
Bombastus	Anionic hydrophilic cream SR...	302200	4
Fagron	Anionic hydrophilic cream SR...	302382	1
Caelo	Anionic hydrophilic cream SR...	302200	1
Bombastus	Anionic hydrophilic cream DAB	302453	4
Bombastus	Anionic hydrophilic cream DAB	302454	4
Caelo	Anionic hydrophilic cream DAB	302453	1
Caelo	Anionic hydrophilic cream DAB	302456	1
Bombastus	Anionic hydrophilic cream DAB	302456	2
Bombastus	Anionic hydrophilic cream DAB	302461	3
Bombastus	Anionic hydrophilic cream DAB	302462	3
Bombastus	Anionic hydrophilic cream DAB	302463	3
Bombastus	Anionic hydrophilic cream DAB	302464	1
Bombastus	Anionic hydrophilic cream DAB	302466	3
Bombastus	Anionic hydrophilic cream DAB	302467	2
Bombastus	Anionic hydrophilic cream DAB	302469	1
Bombastus	Anionic hydrophilic cream SR...	303726	1
Bombastus	Anionic hydrophilic cream DAB	302470	1
Bombastus	Anionic hydrophilic cream SR...	303728	3
Bombastus	Anionic hydrophilic cream SR...	303729	1
Bombastus	Anionic hydrophilic cream SR...	303731	4
Bombastus	Anionic hydrophilic cream SR...	303732	4
Bombastus	Anionic hydrophilic cream SR...	303734	3
Bombastus	Anionic hydrophilic cream SR...	303733	2
Bombastus	Anionic hydrophilic cream SR...	303735	6
Alte Apotheke	Anionic hydrophilic cream SR...	303735	1
Henry Lamotte / Noweda	Anionic hydrophilic cream SR...	303735	2
Bombastus	Anionic hydrophilic cream SR...	303736	3
Bombastus	Anionic hydrophilic cream DAB	304089	6
Bombastus	Anionic hydrophilic cream DAB	304093	1
Bombastus	Anionic hydrophilic cream DAB	304094	1
Caelo	Anionic hydrophilic cream DAB	304096	1
Bombastus	Anionic hydrophilic cream DAB	304098	3
Bombastus	Anionic hydrophilic cream DAB	304097	1
Bombastus	Anionic hydrophilic cream DAB	304099	3
Bombastus	Anionic hydrophilic cream DAB	304100	1
Bombastus	Anionic hydrophilic cream DAB	304101	2
Bombastus	Anionic hydrophilic cream DAB	304102	2
Bombastus	Anionic hydrophilic cream SR...	305064	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	304105	1
Phönix	Anionic hydrophilic cream SR...	305064	1
Bombastus	Anionic hydrophilic cream SR...	305065	2
Bombastus	Anionic hydrophilic cream SR...	305066	2
Bombastus	Anionic hydrophilic cream DAB	305067	1
Phönix	Anionic hydrophilic cream SR...	305066	1
Bombastus	Anionic hydrophilic cream SR...	305067	1
Bombastus	Anionic hydrophilic cream SR...	305068	8
Bombastus	Anionic hydrophilic cream SR...	305069	8
Caelo	Anionic hydrophilic cream SR...	305069	1
Phönix	Anionic hydrophilic cream SR...	305069	1
Frago/Gehe	Anionic hydrophilic cream DAB	305368	5
Bombastus	Anionic hydrophilic cream DAB	305371	1
Euro OTC	Anionic hydrophilic cream DAB	305372	1
Bombastus	Anionic hydrophilic cream DAB	305369	1
Caelo	Anionic hydrophilic cream DAB	305372	1
Bombastus	Anionic hydrophilic cream DAB	305373	2
Bombastus	Anionic hydrophilic cream SR...	305375	1
Bombastus	Anionic hydrophilic cream DAB	305378	3
Bombastus	Anionic hydrophilic cream DAB	305379	2
Bombastus	Anionic hydrophilic cream DAB	305384	1
Bombastus	Anionic hydrophilic cream DAB	305385	2
Bombastus	Anionic hydrophilic cream SR...	306435	7
Bombastus	Anionic hydrophilic cream SR...	306434	1
Bombastus	Anionic hydrophilic cream SR...	306436	1
Bombastus	Anionic hydrophilic cream SR...	306437	1
Bombastus	Anionic hydrophilic cream SR...	306439	1
Bombastus	Anionic hydrophilic cream SR...	306440	4
Bombastuts/Noweda	Anionic hydrophilic cream DAB	307023	2
Bombastus	Anionic hydrophilic cream DAB	307025	2
Bombastus	Anionic hydrophilic cream DAB	307027	2
Bombastus	Anionic hydrophilic cream DAB	307029	4
Bombastus	Anionic hydrophilic cream DAB	307030	2
Bombastus	Anionic hydrophilic cream DAB	307035	1
Bombastus	Anionic hydrophilic cream DAB	307032	1
Bombastus	Anionic hydrophilic cream DAB	307036	5
Bombastus	Anionic hydrophilic cream DAB	307037	2
Bombastus	Anionic hydrophilic cream SR...	307322	1
Bombastus	Anionic hydrophilic cream SR...	307328	1
Bombastus	Anionic hydrophilic cream SR...	307324	1
Bombastus	Anionic hydrophilic cream SR...	307329	4
Bombastus	Anionic hydrophilic cream SR...	307331	2
Bombastus	Anionic hydrophilic cream SR...	307333	3
Bombastus	Anionic hydrophilic cream SR...	307334	2
Fagron	Anionic hydrophilic cream SR...	307334	1
Bombastus	Anionic hydrophilic cream SR...	308214	1
Bombastus	Anionic hydrophilic cream SR...	308217	2
Bombastus	Anionic hydrophilic cream SR...	308218	1
Bombastus	Anionic hydrophilic cream SR...	308220	1
Bombastus	Anionic hydrophilic cream SR...	308221	2
Bombastuts/Gehe	Anionic hydrophilic cream DAB	308397	4
Bombastus	Anionic hydrophilic cream DAB	308401	1
Bombastus	Anionic hydrophilic cream DAB	308402	1
Bombstus/San	Anionic hydrophilic cream DAB	308403	1
Bombastus	Anionic hydrophilic cream SR...	309912	2
apomix	Anionic hydrophilic cream SR...	30A160816	1
Bombastus	Anionic hydrophilic cream DAB	310217	2
Bombastus	Anionic hydrophilic cream DAB	310218	1
Bombastus	Anionic hydrophilic cream DAB	310221	6
Euro OTC	Anionic hydrophilic cream DAB	310266	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	310272	1
Bombastus	Anionic hydrophilic cream SR...	311053	1
Bombastus	Anionic hydrophilic cream SR...	311050	1
Bombastus	Anionic hydrophilic cream DAB	311054	1
Bombastus	Anionic hydrophilic cream SR...	311054	2
Bombastus	Anionic hydrophilic cream SR...	311055	1
Caelo	Anionic hydrophilic cream SR...	311055	1
Bombastus	Anionic hydrophilic cream DAB	311118	4
Bombastus	Anionic hydrophilic cream DAB	311121	5
Bombastus	Anionic hydrophilic cream DAB	311122	6
Euro OTC	Anionic hydrophilic cream DAB	311124	1
Bombastus	Anionic hydrophilic cream DAB	311134	3
Bombastus	Anionic hydrophilic cream DAB	311124	1
AMH Niemann / G	Anionic hydrophilic cream SR...	31A230914	1
Noweda/Apomix	Anionic hydrophilic cream SR...	31A14014	1
GEHE	Anionic hydrophilic cream DAB	311134	1
apomix/Gehe	Anionic hydrophilic cream SR...	31A230914	2
Alliance Healthcare/Apo...	Anionic hydrophilic cream SR...	31A230914	1
Caelo	Anionic hydrophilic cream SR...	31A230914	1
apomix/AHD	Anionic hydrophilic cream SR...	32A150915	1
Bombastus	Anionic hydrophilic cream DAB	320461	1
Apomix	Anionic hydrophilic cream SR...	33A060916	1
Apomix/Noweda	Anionic hydrophilic cream SR...	33A141014	1
Apomix/Gehe	Anionic hydrophilic cream SR...	33A171017	1
Caelo	Anionic hydrophilic cream SR...	33A171017	1
apomix/Phönix	Anionic hydrophilic cream SR...	34A130916	1
AMH Niemann / G	Anionic hydrophilic cream SR...	33A230915	1
Allianc. Health. / Apomix	Anionic hydrophilic cream SR...	34A130916	1
Allianc. Healthc. /Apomix	Anionic hydrophilic cream SR...	34A131015	2
Alliance Healthcare/Apo...	Anionic hydrophilic cream SR...	35A031114	1
Anzag/apomix/07.01.15	Anionic hydrophilic cream SR...	35A031114	1
Apomix/Noweda	Anionic hydrophilic cream SR...	35A140916	1
apomix / Noweda	Anionic hydrophilic cream SR...	36A051217	1
AMH Niemann/G	Anionic hydrophilic cream SR...	36A031115	1
Apomix/Gehe	Anionic hydrophilic cream SR...	36A051217	2
Phönix	Anionic hydrophilic cream SR...	36a140916	1
apomix/Phönix	Anionic hydrophilic cream SR...	37A171115	1
Allianc. Health. /Apomix	Anionic hydrophilic cream SR...	38A111016	2
Apomix/Noweda	Anionic hydrophilic cream SR...	37A171115	1
Phönix	Anionic hydrophilic cream SR...	38a111016	1
PKH Apomix/Anzag	Anionic hydrophilic cream SR...	38A171115	1
Noweda	Anionic hydrophilic cream SR...	38A171115	1
Apomix	Anionic hydrophilic cream SR...	39A011214	2
Rabenapotheke	Anionic hydrophilic cream SR...	39A011214	1
Allianc. Health. /Apomix	Anionic hydrophilic cream SR...	39A231115	1
Bombastus	Anionic hydrophilic cream DAB	403099	1
PKH Halle / Phoenix	Anionic hydrophilic cream SR...	40A091214	1
Apomix	Anionic hydrophilic cream SR...	40A151116	1
apomix/Phönix	Anionic hydrophilic cream SR...	40A151116	1
apomix/Gehe	Anionic hydrophilic cream SR...	41A101214	1
AMH Niemann / G	Anionic hydrophilic cream SR...	41A101214	1
Alliance Healthcare/Apo...	Anionic hydrophilic cream SR...	41A101214	1
Allianc. Health. /Apomix	Anionic hydrophilic cream SR...	41A151215	1
Allianc. Health. / Apomix	Anionic hydrophilic cream SR...	41A151215a	1
Apomix/Sanacorp	Anionic hydrophilic cream SR...	41A281116	1
Apomix/Noweda	Anionic hydrophilic cream SR...	42A171214	1
Apomix	Anionic hydrophilic cream SR...	42A171214	2
Gehe 10.03.15	Anionic hydrophilic cream SR...	42A171214	1
AMH Niemann/E	Anionic hydrophilic cream SR...	45A141113	1
Gehe/Apomix	Anionic hydrophilic cream SR...	48A031213	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream wi...	56/AKD15	1
Bombastus	Anionic hydrophilic cream DAB	60080601	1
Caelo	Anionic hydrophilic cream wi...	64/AKD16	1
Caelo	Anionic hydrophilic cream DAB	9/AKD15	1
Bombastus	Anionic hydrophilic cream SR...	91770	1
Bombastus	Anionic hydrophilic cream SR...	98245	1
Caelo	Anionic hydrophilic cream DAB		8
Caelo	Anionic hydrophilic cream wit...		1
Caelo	Anionic hydrophilic cream wi...		1
Caelo	Anionic hydrophilic cream DAB	Mh090716Lo1	1
Bombastus	Anionic hydrophilic cream SR...	244379	1
Bombastus	Anionic hydrophilic cream SR...	229317	1
Bombastus	Anionic hydrophilic cream SR...	261009	1
Bombastus	Anionic hydrophilic cream SR...	258517	1
Bombastus	Anionic hydrophilic cream SR...	261013	2
Bombastus	Anionic hydrophilic cream SR...	264751	3
Bombastus	Anionic hydrophilic cream SR...	266051	2
Gehe	Anionic hydrophilic cream SR...	267640	1
Bombastus	Anionic hydrophilic cream SR...	267642	2
Bombastus	Anionic hydrophilic cream SR...	269225	1
Bombastus	Anionic hydrophilic cream SR...	269227	2
Gehe	Anionic hydrophilic cream SR...	269226	1
Bombastus	Anionic hydrophilic cream SR...	269231	1
Bombastus	Anionic hydrophilic cream SR...	269234	1
Bombastus	Anionic hydrophilic cream SR...	269235	1
Bombastus	Anionic hydrophilic cream SR...	269236	2
Bombastus	Anionic hydrophilic cream SR...	269317	1
Bombastus	Anionic hydrophilic cream SR...	271701	2
Defektur	Anionic hydrophilic cream SR...	271711	1
Gehe	Anionic hydrophilic cream SR...	271838	1
Noweda	Anionic hydrophilic cream SR...	271711	1
Bombastus	Anionic hydrophilic cream SR...	271837	2
Phönix	Anionic hydrophilic cream SR...	272515	3
Bombastus	Anionic hydrophilic cream SR...	272516	1
Bombastus	Anionic hydrophilic cream SR...	273694	2
Bombastus	Anionic hydrophilic cream SR...	273693	1
Gehe	Anionic hydrophilic cream SR...	273694	1
Bombastus	Anionic hydrophilic cream SR...	273695	2
Noweda	Anionic hydrophilic cream SR...	273695	1
Sanacorp	Anionic hydrophilic cream SR...	274376	1
Caelo	Anionic hydrophilic cream SR...	274377	2
Anzag	Anionic hydrophilic cream SR...	274377	1
Bombastus	Anionic hydrophilic cream SR...	274377	1
Bombastus	Anionic hydrophilic cream SR...	274378	3
Bombastus	Anionic hydrophilic cream SR...	274379	1
AMH Niemann/E	Anionic hydrophilic cream SR...	1A080113	3
Apomix	Anionic hydrophilic cream wit...	2B250213	1
Apomix / Phönix	Anionic hydrophilic cream SR...	3A280114	1
Apomix	Anionic hydrophilic cream SR...	3A150113	1
AMH Niemann GmbH & Co K...	Anionic hydrophilic cream SR...	7A270213	2
Apomix	Anionic hydrophilic cream SR...	5A200212	1
Caelo	Anionic hydrophilic cream SR...	7A2702213	1
Apomix	Anionic hydrophilic cream wit...	6B220513	1
AMH Niemann / Phoenix	Anionic hydrophilic cream wit...	7B180613	1
Caelo	Anionic hydrophilic cream SR...	8061203	1
Apomix	Anionic hydrophilic cream wit...	7B180613	1
Bombastus	Anionic hydrophilic cream SR...	80616-09	1
apomix/GEHE/17.04.2013	Anionic hydrophilic cream SR...	8A270213	1
apomix.de	Anionic hydrophilic cream SR...	8A270213	1
apomix	Anionic hydrophilic cream SR...	8A290311	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream SR...	9011304	1
AMH Niemann/G	Anionic hydrophlic cream wit...	10B160712	1
Caelo	Anionic hydrophilic cream SR...	11101301	1
Pönix/AMH	Anionic hydrophilic cream SR...	11A070313	1
Bombastus	Anionic hydrophilic cream SR...	12121201	1
Caelo	Anionic hydrophilic cream SR...	13072401	1
Phönix Gotha	Anionic hydrophilic cream SR...	122131	3
Apomix	Anionic hydrophlic cream wit...	12B040912	1
Caelo	Anionic hydrophilic cream SR...	13270103	2
Caelo	Anionic hydrophilic cream SR...	13128701	1
apomix/E	Anionic hydrophilic cream SR...	13A190313	1
cae	Anionic hydrophilic cream SR...	14182101	1
	Anionic hydrophilic cream SR...	14182102	1
Caelo	Anionic hydrophilic cream SR...	14182104	1
PKH GmbH/ Phönix	Anionic hydrophilic cream SR...	14A030413	1
Apomix/Alliance Healthc...	Anionic hydrophilic cream SR...	15A040413	3
AMH Niemann/K	Anionic hydrophlic cream wit...	14B161213	1
Caelo	Anionic hydrophilic cream SR...	15011413	1
Caelo	Anionic hydrophilic cream SR...	15069404	1
Caelo	Anionic hydrophilic cream SR...	15111306	1
Bombastus	Anionic hydrophilic cream SR...	15121201	1
Apomix	Anionic hydrophilic cream SR...	15A040413	1
AMH Niemann GmbH&Co.KG	Anionic hydrophilic cream SR...	15A040413	1
apomix/GEHE/10.07.2013	Anionic hydrophilic cream SR...	15A040413	1
PKH Halle / Phoenix	Anionic hydrophilic cream SR...	16A240413	2
Apomix	Anionic hydrophilic cream SR...	16A060611	1
apomix	Anionic hydrophilic cream SR...	16A240413	3
Apomix	Anionic hydrophilic cream SR...	17A240413	2
AMH Niemann/Anzag	Anionic hydrophilic cream SR...	17A240413	1
AMH Niemann/K	Anionic hydrophlic cream wit...	17B271212	1
Bombastus	Anionic hydrophilic cream SR...	18A070513	1
Caelo	Anionic hydrophilic cream SR...	18121202	1
Phoenix	Anionic hydrophilic cream SR...	18A260514	1
AMH Niemann/G	Anionic hydrophilic cream SR...	18A180612	1
Apomix	Anionic hydrophilic cream SR...	18A180612	1
apomix	Anionic hydrophilic cream SR...	18A180612	1
Apomix	Anionic hydrophilic cream SR...	19A070513	1
Caelo	Anionic hydrophilic cream SR...	20061205	1
Caelo	Anionic hydrophilic cream SR...	2013111213	1
Eigene Herstellung	Anionic hydrophilic cream SR...	2014030214	1
Caelo	Anionic hydrophilic cream SR...	21111202	1
Caelo	Anionic hydrophilic cream SR...	22061205	1
Caelo	Anionic hydrophilic cream SR...	20A080513	1
Caelo	Anionic hydrophilic cream SR...	23041307	1
Phönix Gotha	Anionic hydrophilic cream SR...	241111	1
Apomix	Anionic hydrophilic cream SR...	25A050613	1
Caelo	Anionic hydrophilic cream SR...	25021413	1
AMH Niemann/G	Anionic hydrophilic cream SR...	24A070812	1
ApoMix	Anionic hydrophilic cream SR...	25A210812	1
Apomix	Anionic hydrophilic cream SR...	26A220812	1
Caelo	Anionic hydrophilic cream SR...	268979	1
Gehe	Anionic hydrophilic cream SR...	272509	1
Bombastus	Anionic hydrophilic cream SR...	272510	1
Bombastus	Anionic hydrophilic cream SR...	271699	1
Bombastus	Anionic hydrophilic cream SR...	272511	2
Bombastus	Anionic hydrophilic cream SR...	272514	1
Gehe	Anionic hydrophilic cream SR...	273692	1
Bombastus	Anionic hydrophilic cream SR...	274373	1
Gehe	Anionic hydrophilic cream SR...	274375	1
Bombastus	Anionic hydrophilic cream SR...	274381	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fiebi 23.11.12 ek:4,70	Anionic hydrophilic cream SR...	274504	2
Bombastus	Anionic hydrophilic cream SR...	274504	2
Bombastus	Anionic hydrophilic cream SR...	274505	1
Bombastus	Anionic hydrophilic cream SR...	275926	2
Bombastus	Anionic hydrophilic cream SR...	275927	3
Bombastus	Anionic hydrophilic cream SR...	276073	5
Bombastus	Anionic hydrophilic cream SR...	276074	8
Noweda	Anionic hydrophilic cream SR...	276074	1
Euro OTC	Anionic hydrophilic cream SR...	276074	1
Bombastus	Anionic hydrophilic cream SR...	276075	3
Caelo	Anionic hydrophilic cream SR...	276077	1
Bombastus	Anionic hydrophilic cream SR...	277560	2
Bombastus	Anionic hydrophilic cream SR...	277563	2
Gehe	Anionic hydrophilic cream SR...	277562	1
gehe	Anionic hydrophilic cream SR...	277563	1
Bombastus	Anionic hydrophilic cream SR...	277565	3
Caelo	Anionic hydrophilic cream SR...	277565	1
Gehe	Anionic hydrophilic cream SR...	277565	1
Bombastus	Anionic hydrophilic cream SR...	277566	3
Bombastus	Anionic hydrophilic cream SR...	277569	3
Bombastus	Anionic hydrophilic cream SR...	277567	2
Bombastus	Anionic hydrophilic cream SR...	277568	3
Anzag ek: 5,17EUR 08.04.13	Anionic hydrophilic cream SR...	277636	3
Bombastus	Anionic hydrophilic cream SR...	277636	3
Caelo	Anionic hydrophilic cream SR...	277636	1
Bombastus	Anionic hydrophilic cream SR...	277637	5
	Anionic hydrophilic cream SR...	277637	2
Caelo	Anionic hydrophilic cream SR...	277637	1
Bombastus	Anionic hydrophilic cream SR...	277637	1
Sanacorp	Anionic hydrophilic cream SR...	277637	2
Bombastus	Anionic hydrophilic cream SR...	278026	2
Bombastus	Anionic hydrophilic cream SR...	279195	2
Bombastus	Anionic hydrophilic cream SR...	279198	2
Bombastus	Anionic hydrophilic cream SR...	279199	3
Bombastus	Anionic hydrophilic cream SR...	279200	2
Bombastus	Anionic hydrophilic cream SR...	279201	2
Bombastus	Anionic hydrophilic cream SR...	279323	3
Bombastus	Anionic hydrophilic cream SR...	279324	3
Bombastus	Anionic hydrophilic cream SR...	279325	10
Bombastus	Anionic hydrophilic cream SR...	279326	10
Noweda	Anionic hydrophilic cream SR...	279326	1
Apomix	Anionic hydrophilic cream SR...	27A180613	1
Apomix	Anionic hydrophilic cream SR...	27A210911	1
Bombastus	Anionic hydrophilic cream SR...	280570	2
Noweda	Anionic hydrophilic cream SR...	280572	1
Bombastus	Anionic hydrophilic cream SR...	280571	1
Bombastus	Anionic hydrophilic cream SR...	280572	1
Bombastus	Anionic hydrophilic cream SR...	280574	4
Gehe	Anionic hydrophilic cream SR...	280574	1
Bombastus	Anionic hydrophilic cream SR...	280575	1
Gehe	Anionic hydrophilic cream SR...	280575	1
Bombastus	Anionic hydrophilic cream SR...	280576	10
Bombastus	Anionic hydrophilic cream SR...	280741	12
Noweda	Anionic hydrophilic cream SR...	280741	2
Bombastus	Anionic hydrophilic cream SR...	281737	3
Caelo	Anionic hydrophilic cream SR...	281738	2
Bombastus	Anionic hydrophilic cream SR...	281739	7
Caelo	Anionic hydrophilic cream SR...	281739	3
Bombastus	Anionic hydrophilic cream SR...	281742	4
Bombastus	Anionic hydrophilic cream SR...	281744	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream SR...	281745	1
Bombastus	Anionic hydrophilic cream SR...	281809	4
Bombastus	Anionic hydrophilic cream SR...	283282	3
Bombastus	Anionic hydrophilic cream SR...	284077	1
Alliance Healthcare/Apo...	Anionic hydrophilic cream SR...	28A250613	2
Pharmachem/Phönix	Anionic hydrophilic cream SR...	290413	1
Sanacorp	Anionic hydrophilic cream SR...	290413	2
Caelo	Anionic hydrophilic cream SR...	290530	1
Apomix	Anionic hydrophilic cream SR...	29A240912	1
Bombastus	Anionic hydrophilic cream SR...	301357	2
All..Healthc./Apomix	Anionic hydrophilic cream SR...	30A020713	1
Gehe/apomix/18.09.13	Anionic hydrophilic cream SR...	30A020713	2
AMH Niemann GmbH & Co. KG	Anionic hydrophilic cream SR...	30A0207123	1
Phoenix (PKH Halle)	Anionic hydrophilic cream SR...	30A021012	1
Bombastus	Anionic hydrophilic cream SR...	312121-3-P	1
Caelo	Anionic hydrophilic cream SR...	31071301	1
Apomix	Anionic hydrophilic cream SR...	31A030713	1
apomix	Anionic hydrophilic cream SR...	31A030713	1
AMH Niemann/E	Anionic hydrophilic cream SR...	31A091012	1
Apomix	Anionic hydrophilic cream SR...	32A160713	1
Apomix	Anionic hydrophilic cream SR...	33A220713	1
Phoenix (PKH Halle)	Anionic hydrophilic cream SR...	33A151012	1
Caelo	Anionic hydrophilic cream SR...	34A060813	2
Apomix/Phoenix	Anionic hydrophilic cream SR...	33A220713	1
Apomix	Anionic hydrophilic cream SR...	34A060813	1
GEHE/apomix/01.03.2012	Anionic hydrophilic cream SR...	33A291111	1
AMH Niemann/K	Anionic hydrophilic cream SR...	36A270813	1
AMH Niemann/G	Anionic hydrophilic cream SR...	38A271112	1
Eigenherstellung	Anionic hydrophilic cream SR...	3812E-07664	1
Apomix	Anionic hydrophilic cream SR...	41A021013	1
apomix	Anionic hydrophilic cream SR...	35A081112	1
AMH Niemann/E	Anionic hydrophilic cream SR...	42A141013	1
Eigenherstellung	Anionic hydrophilic cream SR...	4012Q-07664	1
AMH Niemann/G	Anionic hydrophilic cream SR...	44A051113	1
Apomix	Anionic hydrophilic cream SR...	44A051113	1
Apomix	Anionic hydrophilic cream SR...	4546279	1
Apomix	Anionic hydrophilic cream SR...	48A031213	1
Anzag/apomix/31.01.2014	Anionic hydrophilic cream SR...	51A171213	1
Bombastus	Anionic hydrophilic cream SR...	5012E-03327	1
Bombastus	Anionic hydrophilic cream SR...	60089211	1
Bombastus	Anionic hydrophilic cream SR...	82440	1
Caelo	Base cream, hydrophobic DAC	4004319	1
Pharmachem	Base cream, hydrophobic DAC	110413	7
Caesar & Loretz GmbH Hi...	Base cream, hydrophobic DAC	115421302	1
Caelo	Base cream, hydrophobic DAC	13012015A	1
Fagron	Base cream, hydrophobic DAC	13H13-T06	2
Caelo	Base cream, hydrophobic DAC	13H13-T06-006165	2
Caelo	Base cream, hydrophobic DAC	14000701	19
Caelo	Base cream, hydrophobic DAC	1400701	1
Caelo	Base cream, hydrophobic DAC	14028101	26
Phoenix	Base cream, hydrophobic DAC	14028101	1
Caelo	Base cream, hydrophobic DAC	14028102	10
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	140318	1
Caelo	Base cream, hydrophobic DAC	14061402	5
Caelo	Base cream, hydrophobic DAC	14061501	24
Euro OTC	Base cream, hydrophobic DAC	14061501	1
Phönix	Base cream, hydrophobic DAC	14061501	1
Caelo	Base cream, hydrophobic DAC	14061601	5
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	14061501	3
Pharmachen	Base cream, hydrophobic DAC	140885	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PHARMACHEM / 20	Base cream, hydrophobic DAC	140885	1
gehe	Base cream, hydrophobic DAC	140886	2
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	140886	1
Pharmachem/ Sanacorp	Base cream, hydrophobic DAC	140886	1
Caelo	Base cream, hydrophobic DAC	14120903	4
Pharmachem	Base cream, hydrophobic DAC	141317	2
Pharmachem.Pöfneck/Phoe...	Base cream, hydrophobic DAC	141317	1
Gehe	Base cream, hydrophobic DAC	141317	2
Pharmachem/30	Base cream, hydrophobic DAC	141317	2
Gehe	Base cream, hydrophobic DAC	141318	2
Pharmachem Pösneck/ Now...	Base cream, hydrophobic DAC	141318	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	141318	1
Pharmachem Anzag	Base cream, hydrophobic DAC	141318	1
Caelo	Base cream, hydrophobic DAC	14145202	5
Caelo	Base cream, hydrophobic DAC	14145301	20
Phoenix	Base cream, hydrophobic DAC	14145401	2
AHD	Base cream, hydrophobic DAC	14145401	1
Caelo	Base cream, hydrophobic DAC	14145401	23
Berkel AHK	Base cream, hydrophobic DAC	14145401	1
Pharmachem GmbH&Co.KG	Base cream, hydrophobic DAC	141489	1
Kehr / Pharmachem	Base cream, hydrophobic DAC	141489	1
Pharmachem	Base cream, hydrophobic DAC	142154	4
Caelo	Base cream, hydrophobic DAC	14207002	3
Pharmachem/ Noweda	Base cream, hydrophobic DAC	142154	1
Noweda/Pharmachem	Base cream, hydrophobic DAC	142154	1
Pharmachem	Base cream, hydrophobic DAC	142155	2
Pharmachem/Phoenix	Base cream, hydrophobic DAC	142154	2
Pharamachem/Alliance He...	Base cream, hydrophobic DAC	142154	1
Gehe	Base cream, hydrophobic DAC	142155	1
Caelo	Base cream, hydrophobic DAC	14215701	11
Pharmachem/Alliance Hea...	Base cream, hydrophobic DAC	142155	1
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	143093	2
Pharmachem	Base cream, hydrophobic DAC	143093	6
Pharmachem Pöfneck/Sana...	Base cream, hydrophobic DAC	143093	2
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	143094	1
Pharmachem	Base cream, hydrophobic DAC	143094	6
Pharmachen GmbH/Gehe	Base cream, hydrophobic DAC	143094	3
Gehe	Base cream, hydrophobic DAC	143119	1
Pharmachen	Base cream, hydrophobic DAC	143094	1
Caelo	Base cream, hydrophobic DAC	14335501	14
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	14335501	2
Caelo	Base cream, hydrophobic DAC	14335502	9
AH 24.02.15	Base cream, hydrophobic DAC	14335502	1
Caelo	Base cream, hydrophobic DAC	14335601	6
Caelo	Base cream, hydrophobic DAC	14379902	4
Fagron	Base cream, hydrophobic DAC	14A06-T07-009033	3
Fagron	Base cream, hydrophobic DAC	14A06-T08-009037	10
Fagron	Base cream, hydrophobic DAC	14A06T08009038	1
Fagron	Base cream, hydrophobic DAC	14L18-T21-018127	1
Fagron	Base cream, hydrophobic DAC	14L18-T21-018124	2
Caelo	Base cream, hydrophobic DAC	15001001	8
Caelo	Base cream, hydrophobic DAC	15001002	10
Caelo	Base cream, hydrophobic DAC	15001003	2
Pharmachem	Base cream, hydrophobic DAC	150404	3
Pharmachem/AHD	Base cream, hydrophobic DAC	150405	1
Caelo	Base cream, hydrophobic DAC	15030702	1
Caelo	Base cream, hydrophobic DAC	150405	1
Pharmachen Pöfneck	Base cream, hydrophobic DAC	150405	2
Caelo	Base cream, hydrophobic DAC	15044702	6
Caelo	Base cream, hydrophobic DAC	15046401	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream, hydrophobic DAC	15046402	6
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	15048701	2
Caelo	Base cream, hydrophobic DAC	15048701	13
Caelo	Base cream, hydrophobic DAC	15048702	3
Caelo	Base cream, hydrophobic DAC	15111401	7
Caelo	Base cream, hydrophobic DAC	15111402	15
Pharmachem	Base cream, hydrophobic DAC	151257	4
Pharmachem Pößneck/Sana. ...	Base cream, hydrophobic DAC	151257	1
Anzag	Base cream, hydrophobic DAC	151257	1
Pharmachem Anzag	Base cream, hydrophobic DAC	151257	1
Pharmachen GmbH/Gehe	Base cream, hydrophobic DAC	151258	2
Pößneck Pharma / Gehe	Base cream, hydrophobic DAC	151258	1
Pharmachen	Base cream, hydrophobic DAC	151258	4
Pharmachem/Gehe	Base cream, hydrophobic DAC	151259	5
Caelo	Base cream, hydrophobic DAC	151259	3
Pharmachem	Base cream, hydrophobic DAC	151259	3
Caelo	Base cream, hydrophobic DAC	15132502	5
Pharmachem/ Noweda	Base cream, hydrophobic DAC	151259	1
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	151259	2
Caelo	Base cream, hydrophobic DAC	1513502	1
Pharmachen GmbH	Base cream, hydrophobic DAC	152209	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	152209	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	152209	3
Pharmachem/ Noweda	Base cream, hydrophobic DAC	152209	1
Pharmachen	Base cream, hydrophobic DAC	152209	1
Noweda	Base cream, hydrophobic DAC	152209	1
Pharmachem	Base cream, hydrophobic DAC	152209	3
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	15221	1
AHD/ Pharmachem	Base cream, hydrophobic DAC	152210	2
Pharmachem Pößneck/Sana. ...	Base cream, hydrophobic DAC	152210	1
Pharmachem	Base cream, hydrophobic DAC	152210	2
Pharmachem	Base cream, hydrophobic DAC	152211	3
Pharmachem/Gehe	Base cream, hydrophobic DAC	152211	1
Caelo	Base cream, hydrophobic DAC	15221202	17
Caelo	Base cream, hydrophobic DAC	15221301	19
Noweda	Base cream, hydrophobic DAC	15221301	4
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	15221401	1
Pharmachem	Base cream, hydrophobic DAC	153625	6
Pharmachem/Gehe	Base cream, hydrophobic DAC	153625	1
Pharmachem Pößneck/Sana. ...	Base cream, hydrophobic DAC	153626	1
Pharmachem /Phönix	Base cream, hydrophobic DAC	153626	1
Pharmachem Pößneck/Sana. ...	Base cream, hydrophobic DAC	153625	1
Pharmachem/ Anzag	Base cream, hydrophobic DAC	153627	1
Pharmachem/Kehr	Base cream, hydrophobic DAC	153627	1
Pharmachem	Base cream, hydrophobic DAC	153627	4
Caelo	Base cream, hydrophobic DAC	15421201	24
WE:06.04.2016 Phönix ...	Base cream, hydrophobic DAC	15421201	1
Cael/Sanacorp	Base cream, hydrophobic DAC	15421302	1
Caelo	Base cream, hydrophobic DAC	15421302	8
Caelo	Base cream, hydrophobic DAC	15421401	13
Fagron	Base cream, hydrophobic DAC	15J13-T01-026887	4
Fagron	Base cream, hydrophobic DAC	15J13-T01-026888	1
pharmachem pößneck	Base cream, hydrophobic DAC	160130	1
Caelo	Base cream, hydrophobic DAC	160130	3
Pharmachen	Base cream, hydrophobic DAC	160130	1
Ebert u. Jacobi 16.0...	Base cream, hydrophobic DAC	160130	1
Pharmachem/Anzag	Base cream, hydrophobic DAC	160136	1
Pharmachem	Base cream, hydrophobic DAC	160130	6
Pharmachem	Base cream, hydrophobic DAC	160136	5
Pharmachem/ Noweda	Base cream, hydrophobic DAC	160136	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream, hydrophobic DAC	160136	1
Gehe-Pharmachem	Base cream, hydrophobic DAC	160136	1
Pharmachem Pöckneck/Sana...	Base cream, hydrophobic DAC	160137	2
AHD /Pharmachem Pöckneck	Base cream, hydrophobic DAC	160137	3
Pharmachem/Alliance Hea...	Base cream, hydrophobic DAC	160138	1
Caelo	Base cream, hydrophobic DAC	160575101	1
Caelo	Base cream, hydrophobic DAC	16074801	6
Caelo	Base cream, hydrophobic DAC	16074802	16
Caelo	Base cream, hydrophobic DAC	16074901	14
Caelo	Base cream, hydrophobic DAC	16074803	5
Noweda	Base cream, hydrophobic DAC	16074901	2
Caelo	Base cream, hydrophobic DAC	16075101	8
Phoenix	Base cream, hydrophobic DAC	16075101	2
Caelo	Base cream, hydrophobic DAC	16075201	12
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	160814	1
Bombastus	Base cream, hydrophobic DAC	160814	1
Pharmachem GmbH & Co. KG	Base cream, hydrophobic DAC	160825	1
Pharmachem Pöckneck	Base cream, hydrophobic DAC	160825	2
Pharmachem	Base cream, hydrophobic DAC	160825	4
Pharmachem GmbH&Co.KG	Base cream, hydrophobic DAC	160825	3
Pharmachem/ Gehe	Base cream, hydrophobic DAC	160826	1
Pharmachem Pöckneck/Sana...	Base cream, hydrophobic DAC	160826	2
Caelo	Base cream, hydrophobic DAC	1615501	1
Caelo	Base cream, hydrophobic DAC	16165403	17
Caelo	Base cream, hydrophobic DAC	16165501	22
Caelo	Base cream, hydrophobic DAC	16165502	4
Caelo	Base cream, hydrophobic DAC	16165701	9
Alte Apotheke	Base cream, hydrophobic DAC	16165701	1
Caelo	Base cream, hydrophobic DAC	16174201	12
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	16174201	2
Pharmachem	Base cream, hydrophobic DAC	162167	3
Pharmachem	Base cream, hydrophobic DAC	162169	3
Pharmachem/ Noweda	Base cream, hydrophobic DAC	162167	1
Pharmachem Anzag	Base cream, hydrophobic DAC	162168	1
Pharmachem	Base cream, hydrophobic DAC	162168	3
Pharmchem/Sanacorp	Base cream, hydrophobic DAC	162169	1
Pharmachen/gehe	Base cream, hydrophobic DAC	162170	1
Pharmachem/ Noweda	Base cream, hydrophobic DAC	162169	1
Caelo	Base cream, hydrophobic DAC	16293601	4
Caelo	Base cream, hydrophobic DAC	16293701	12
	Base cream, hydrophobic DAC	16293901	1
Euro OTC	Base cream, hydrophobic DAC	16293901	2
Caelo	Base cream, hydrophobic DAC	16293902	6
Pharmachem Anzag	Base cream, hydrophobic DAC	163238	1
AHD/ Pharmachem	Base cream, hydrophobic DAC	163238	2
Caelo	Base cream, hydrophobic DAC	16330201	4
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	16330203	1
Caelo	Base cream, hydrophobic DAC	16330203	8
Fagron	Base cream, hydrophobic DAC	16E03-T02-032868	5
Caelo	Base cream, hydrophobic DAC	170008	2
AHD/Pharmachem.	Base cream, hydrophobic DAC	170008	1
Pharmachem	Base cream, hydrophobic DAC	170009	2
Caelo	Base cream, hydrophobic DAC	17073201	2
Caelo	Base cream, hydrophobic DAC	17073103	2
pharmachem pöckneck	Base cream, hydrophobic DAC	171127	3
Caelo	Base cream, hydrophobic DAC	171128	2
AHD /Pharmachem Pöckneck	Base cream, hydrophobic DAC	171128	2
Pharmachem Anzag	Base cream, hydrophobic DAC	171128	1
Caelo	Base cream, hydrophobic DAC	17112903	4
Caelo	Base cream, hydrophobic DAC	17150501	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Pharmachen GmbH	Base cream, hydrophobic DAC	171817	1
Pharmachem	Base cream, hydrophobic DAC	171815	2
Sanacorp/Pharmachem	Base cream, hydrophobic DAC	171816	1
Pharmachem GmbH & Co. K...	Base cream, hydrophobic DAC	17185	1
Pharmachem	Base cream, hydrophobic DAC	172194	3
Caelo	Base cream, hydrophobic DAC	17255302	1
Fagron	Base cream, hydrophobic DAC	17C06-T08-042287	2
Caelo	Base cream, hydrophobic DAC	18051101	3
Pharmachem/Gehe	Base cream, hydrophobic DAC	180517	2
Caelo	Base cream, hydrophobic DAC	230115/2	1
gehe	Base cream, hydrophobic DAC	231213	3
Kehr/ Pharmachem	Base cream, hydrophobic DAC	231213	2
Caelo	Base cream, hydrophobic DAC	282163	1
Bombastus	Base cream, hydrophobic DAC	286514	2
Caelo	Base cream, hydrophobic DAC	286516	1
Bombastus	Base cream, hydrophobic DAC	286516	6
Bombastus	Base cream, hydrophobic DAC	288499	4
Bombastus	Base cream, hydrophobic DAC	289718	1
Gehe/apomix/17.12.2014	Base cream, hydrophobic DAC	289719	1
Gehe	Base cream, hydrophobic DAC	289721	1
Bombastus	Base cream, hydrophobic DAC	289720	2
Bombastus	Base cream, hydrophobic DAC	291590	4
Bombastus	Base cream, hydrophobic DAC	291593	3
Bombastus	Base cream, hydrophobic DAC	293168	4
Bombastus	Base cream, hydrophobic DAC	293163	5
Caelo	Base cream, hydrophobic DAC	293163	1
Caelo	Base cream, hydrophobic DAC	294686	3
Bombastus	Base cream, hydrophobic DAC	297043	2
Bombastus	Base cream, hydrophobic DAC	297045	1
Bombastus	Base cream, hydrophobic DAC	298350	1
Bombastus	Base cream, hydrophobic DAC	300352	1
Euro OTC	Base cream, hydrophobic DAC	300415	1
Bombastus	Base cream, hydrophobic DAC	300417	4
Bombastus	Base cream, hydrophobic DAC	301455	2
Bombastus	Base cream, hydrophobic DAC	301457	6
Bombastus	Base cream, hydrophobic DAC	302131	1
Caelo	Base cream, hydrophobic DAC	303376	3
Bombastus	Base cream, hydrophobic DAC	303376	5
Bombastus	Base cream, hydrophobic DAC	305083	5
Bombastus	Base cream, hydrophobic DAC	305080	1
Bombastus	Base cream, hydrophobic DAC	306395	2
Bombastus	Base cream, hydrophobic DAC	306397	4
Bombastus	Base cream, hydrophobic DAC	307771	5
Bombastus	Base cream, hydrophobic DAC	307766	2
Caelo	Base cream, hydrophobic DAC	308507	1
Bombastus	Base cream, hydrophobic DAC	308508	2
Bombastus	Base cream, hydrophobic DAC	310827	4
Pharmachem/Gehe	Base cream, hydrophobic DAC	21013	1
Bombastus	Base cream, hydrophobic DAC	264854	1
Gehe	Base cream, hydrophobic DAC	21013	1
Pharmachen GmbH/Gehe	Base cream, hydrophobic DAC	21013	1
Pharmachem GmbH&Co.KG	Base cream, hydrophobic DAC	21013	1
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	21013	1
Pharmachem	Base cream, hydrophobic DAC	21013	1
Gehe	Base cream, hydrophobic DAC	90512	1
PharmaChem	Base cream, hydrophobic DAC	90512	1
Sanacorp/ Pharmachem	Base cream, hydrophobic DAC	90512	3
Caelo	Base cream, hydrophobic DAC	11038802	2
Caelo	Base cream, hydrophobic DAC	11038904	1
Noweda	Base cream, hydrophobic DAC	110413	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Pharmachem / Anzag	Base cream, hydrophobic DAC	110413	2
Caelo	Base cream, hydrophobic DAC	11222003	1
Caesar & Loretz GmbH, 4. ...	Base cream, hydrophobic DAC	11289203	1
Caelo	Base cream, hydrophobic DAC	11289203	1
Caelo	Base cream, hydrophobic DAC	113401301	1
Caelo	Base cream, hydrophobic DAC	11289204	1
Caelo	Base cream, hydrophobic DAC	11341301	2
Noweda	Base cream, hydrophobic DAC	11341401	1
Caelo	Base cream, hydrophobic DAC	11341402	2
Fagron	Base cream, hydrophobic DAC	11L12-T05	1
Taoasis	Base cream, hydrophobic DAC	120214-2023	1
Noweda	Base cream, hydrophobic DAC	12025002	4
Caelo	Base cream, hydrophobic DAC	12025002	2
Caelo	Base cream, hydrophobic DAC	12025003	2
Taoasis	Base cream, hydrophobic DAC	120601-2137	1
Caelo	Base cream, hydrophobic DAC	12064703	3
Caelo	Base cream, hydrophobic DAC	12064704	5
Pharmachem / Anzag	Base cream, hydrophobic DAC	120712	1
Caelo	Base cream, hydrophobic DAC	12064705	6
Pharmachem GmbH	Base cream, hydrophobic DAC	120712	1
Pharmachem GmbH	Base cream, hydrophobic DAC	120712	2
Taoasis	Base cream, hydrophobic DAC	120719-2199	1
Taoasis	Base cream, hydrophobic DAC	121205-2276	1
Taoasis	Base cream, hydrophobic DAC	121205-2299	1
Noweda	Base cream, hydrophobic DAC	12151601	3
Anzag 30.10.12 ek:12,68EUR	Base cream, hydrophobic DAC	12151601	1
Caelo	Base cream, hydrophobic DAC	12151601	10
Caelo	Base cream, hydrophobic DAC	12151602	7
Caelo	Base cream, hydrophobic DAC	12229302	7
Anzag 07.05.13 ek:12,66EUR	Base cream, hydrophobic DAC	12229302	1
Caelo	Base cream, hydrophobic DAC	12308102	4
Gehe	Base cream, hydrophobic DAC	12308102	1
Caelo	Base cream, hydrophobic DAC	12308201	6
Noweda	Base cream, hydrophobic DAC	12308201	2
Caelo	Base cream, hydrophobic DAC	12308202	2
Fagron	Base cream, hydrophobic DAC	12J24-T01	7
Hirsch	Base cream, hydrophobic DAC	12J24-T01	2
Taoasis	Base cream, hydrophobic DAC	130124-3021	1
Caelo	Base cream, hydrophobic DAC	13019302	6
Caelo	Base cream, hydrophobic DAC	13019303	14
Noweda	Base cream, hydrophobic DAC	13019303	6
Sanacorp	Base cream, hydrophobic DAC	13019303	1
Taoasis	Base cream, hydrophobic DAC	130312-3053	1
Caelo	Base cream, hydrophobic DAC	13019304	3
calo	Base cream, hydrophobic DAC	13072602	3
Caelo	Base cream, hydrophobic DAC	13072602	6
Caelo	Base cream, hydrophobic DAC	13085601	13
Caelo	Base cream, hydrophobic DAC	13085602	8
Caelo	Base cream, hydrophobic DAC	13085701	24
Gehe	Base cream, hydrophobic DAC	131011	1
Kehr	Base cream, hydrophobic DAC	131011	1
Caelo	Base cream, hydrophobic DAC	131211	2
Caelo	Base cream, hydrophobic DAC	13175901	21
Phönix	Base cream, hydrophobic DAC	13175901	3
Gehe	Base cream, hydrophobic DAC	13175901	1
Caelo	Base cream, hydrophobic DAC	13175902	8
Gehe	Base cream, hydrophobic DAC	13175902	1
Caelo	Base cream, hydrophobic DAC	13176003	13
Caelo	Base cream, hydrophobic DAC	13274902	2
Caelo	Base cream, hydrophobic DAC	13275301	12

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Base cream, hydrophobic DAC	13275302	2
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	13275302	2
Caelo	Base cream, hydrophobic DAC	13275303	6
Caelo	Base cream, hydrophobic DAC	13386401	13
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	13386401	2
Gehe	Base cream, hydrophobic DAC	13386402	1
Phönix	Base cream, hydrophobic DAC	13386402	2
Caelo	Base cream, hydrophobic DAC	13B04-T17-001037	2
Caelo	Base cream, hydrophobic DAC	13386403	3
Fagron	Base cream, hydrophobic DAC	13B04-T17-001037	4
Fagron	Base cream, hydrophobic DAC	13E13-T01-003619	6
Fagron	Base cream, hydrophobic DAC	13H13-T06-006165	5
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	14000701	3
Phönix	Base cream, hydrophobic DAC	14000701	1
Phoenix	Base cream, hydrophobic DAC	14061501	2
Sanacorp	Base cream, hydrophobic DAC	140885	1
Caesar&Loretz,Phoenix	Base cream, hydrophobic DAC	14335502	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	230713	2
Caelo	Base cream, hydrophobic DAC	230713	3
Pharmachem GmbH / Anzag	Base cream, hydrophobic DAC	230713	1
Pharmachem	Base cream, hydrophobic DAC	231213	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	240113	2
Pharmachem GmbH	Base cream, hydrophobic DAC	240113	2
Sanacorp/ Pharmachem	Base cream, hydrophobic DAC	261112	2
Bombastus	Base cream, hydrophobic DAC	261112	2
Bombastus	Base cream, hydrophobic DAC	271092	3
Bombastus	Base cream, hydrophobic DAC	276085	1
Phönix Gotha	Base cream, hydrophobic DAC	276087	2
Bombastus	Base cream, hydrophobic DAC	280811	2
Bombastus	Base cream, hydrophobic DAC	279147	2
Bombastus	Base cream, hydrophobic DAC	282162	2
Bombastus	Base cream, hydrophobic DAC	282163	1
Bombastus	Base cream, hydrophobic DAC	284173	1
Bombastus	Base cream, hydrophobic DAC	286515	1
Phönix / InfectoPharm	Base cream, hydrophobic DAC		1
Caelo	Anionic hydrophilic cream SR...	16158901	1
Sanacorp	Anionic hydrophilic cream SR...	293133	2
Bombastus	Anionic hydrophilic cream SR...	295172	1
Caelo	Base cream, hydrophobic DAC	1414501	2
Pharmachem/AHD	Base cream, hydrophobic DAC	151259	2
Pharmachem Pößneck/Phoe...	Base cream, hydrophobic DAC	153627	1
Pharmachem;	Base cream, hydrophobic DAC	160136	1
Caelo	Base cream, hydrophobic DAC	16075003	5
Pharmachem/Phoenix	Base cream, hydrophobic DAC	162168	1
Bombastus	Base cream, hydrophobic DAC	298353	5
Caelo	Base cream, hydrophobic DAC	11153902	1
Caelo	Base cream, hydrophobic DAC	11222002	1
Anzag/Pharmachem	Base cream, hydrophobic DAC	120712	2
Caelo	Base cream, hydrophobic DAC	12292102	1
Gutenberg-Apotheke	Anionic hydrophilic cream SR...	285890	1
Bombastus	Anionic hydrophilic cream SR...	286974	1
Bombastus	Anionic hydrophilic cream SR...	286978	1
Bombastus	Anionic hydrophilic cream SR...	288228	1
VDL;02.12.14;5,32EUR	Anionic hydrophilic cream SR...	288310	1
Caesar & Loretz GmbH	Anionic hydrophilic cream SR...	289329	1
Bombastus	Anionic hydrophilic cream SR...	289324	1
Bombastus	Anionic hydrophilic cream SR...	290526	1
Bombastus	Anionic hydrophilic cream SR...	292973	1
Bombastus	Anionic hydrophilic cream SR...	283281	3
Sanacorp	Anionic hydrophilic cream SR...	283283	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Anionic hydrophilic cream SR...	283284	1
Pharmachem	Skin care ointment W/L SR	20316	1
Apomix	Skin care ointment W/L SR (w...	3A220118	1
Apomix	Skin care ointment W/L SR (w...	4A220118	1
Caelo	Skin care ointment W/L SR	80115	1
Apomix	Skin care ointment W/L SR (w...	8A280218	1
Pharmachem GmbH Pößneck	Skin care ointment W/L SR	81117	1
Pharmachem	Skin care ointment W/L SR	90216	1
Apomix	Skin care ointment W/L SR (w...	9A0103318	1
Pharmachem Pößneck/Noweda	Skin care ointment W/L SR	111016	1
Noweda/Pharmachem	Skin care ointment W/L SR	120416	1
Pharmachem Pößneck	Skin care ointment W/L SR	100516	1
Pharmachem/Noweda	Skin care ointment W/L SR	110416	1
Apomix	Skin care ointment W/L SR	16A250314	2
Apomix	Skin care ointment W/L SR (w...	16A080518	1
Apomix	Skin care ointment W/L SR (w...	15A080518	1
Apomix/ Sanacorp	Skin care ointment W/L SR (w...	15B220216	1
Pharmachem/ Sanacorp	Skin care ointment W/L SR	141015	1
Pharmachem/ Sanacorp	Skin care ointment W/L SR	170117	1
Apomix	Skin care ointment W/L SR	17A250314	1
Pharmachem / Noweda	Skin care ointment W/L SR	170815	1
Pharmachem/ Noweda	Skin care ointment W/L SR	180914	1
Apomix	Skin care ointment W/L SR	21A230414	1
Pharmachem GmbH	Skin care ointment W/L SR	220118	1
Pharmachem/ Noweda	Skin care ointment W/L SR	260814	1
Apomix	Skin care ointment W/L SR	22A230414	1
Noweda/Pharmachem	Skin care ointment W/L SR	270715	1
Euro OTC	Skin care ointment W/L SR	300316	1
Apomix/ Sanacorp	Skin care ointment W/L SR (w...	41B200715	1
apomix	Skin care ointment W/L SR (w...	31B130917	1
pkh gmbh	Skin care ointment W/L SR (w...	49b171014	1
Apomix	Skin care ointment W/L SR	9A260213	1
Apomix	Skin care ointment W/L SR	1A080114	1
Apomix	Skin care ointment W/L SR	2A090114	1
Apomix	Skin care ointment W/L SR	11A190214	1
Apomix	Skin care ointment W/L SR	211112	1
Pharmachem/Gehe	Skin care ointment W/L SR	100712	1
Gehe	Hydrophilic cream, non-ionic...		1
Bombastus	Hydrophilic cream, non-ionic...	229352	1
Apomix	Hydrophilic cream, non-ionic...	1A200116	2
Apomix	Hydrophilic cream, non-ionic...	1B070116	2
apomix	Hydrophilic cream, non-ionic...	3A030216	1
apomix / Anzag	Hydrophilic cream, non-ionic...	3A030216	1
apomix	Hydrophilic cream, non-ionic...	3B020215	1
Apomix	Hydrophilic cream, non-ionic...	5A080316	1
apomix/Phönix	Hydrophilic cream, non-ionic...	5A080316	3
Apomix/Kehr	Hydrophilic cream, non-ionic...	8A120416	1
PKH Halle/Sanacorp	Hydrophilic cream, non-ionic...	8A030414	1
Pharmachem	Hydrophilic cream, non-ionic...	90714	5
Pharmachem GmbH&Co.KG	Hydrophilic cream, non-ionic...	90714	1
apomix	Hydrophilic cream, non-ionic...	9B020516	1
Sanacorp	Anionic hydrophilic cream DAB	104/AKD4	1
Sanacorp	Anionic hydrophilic cream DAB	106/AKD12	1
Apomix / 11	Hydrophilic cream, non-ionic...	10B060514	1
Fagron	Hydrophilic cream, non-ionic...	11A180416	1
Caelo	Anionic hydrophilic cream DAB	12886	1
Caelo	Anionic hydrophilic cream DAB	1259487	1
Phönix Gotha	Hydrophilic cream, non-ionic...	12B211015	1
Apomix	Hydrophilic cream, non-ionic...	12B210514	1
Caelo	Anionic hydrophilic cream DAB	13265501	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	13265602	7
Caelo	Anionic hydrophilic cream DAB	13314201	9
Caelo	Anionic hydrophilic cream DAB	13315001	7
Caelo	Anionic hydrophilic cream DAB	13346201	3
Caelo	Anionic hydrophilic cream DAB	13363301	1
Caelo	Anionic hydrophilic cream DAB	13363201	11
Caelo	Anionic hydrophilic cream DAB	13363601	13
Caelo	Anionic hydrophilic cream DAB	13364001	13
Caelo	Anionic hydrophilic cream DAB	13364101	4
Caelo	Anionic hydrophilic cream DAB	14016701	9
Caelo	Anionic hydrophilic cream DAB	14016901	1
Caelo	Anionic hydrophilic cream DAB	14024901	1
Phönix Gotha	Hydrophilic cream, non-ionic...	140551	4
Caelo	Anionic hydrophilic cream DAB	14061801	9
Caelo	Anionic hydrophilic cream DAB	14062101	7
Großhandel	Anionic hydrophilic cream DAB	14062101	1
Caelo	Anionic hydrophilic cream DAB	14062201	9
Großhandel	Anionic hydrophilic cream DAB	14062301	1
Caelo	Anionic hydrophilic cream DAB	14062301	8
Caelo	Anionic hydrophilic cream DAB	14062701	4
Caelo	Anionic hydrophilic cream DAB	14062601	10
Caelo	Anionic hydrophilic cream DAB	14062901	10
Caelo	Anionic hydrophilic cream DAB	14062801	4
Caelo	Hydrophilic cream, non-ionic...	14063001	2
Caelo	Anionic hydrophilic cream DAB	14068601	1
Caelo	Anionic hydrophilic cream DAB	14118101	15
Caelo	Anionic hydrophilic cream DAB	14118201	11
Caelo	Anionic hydrophilic cream DAB	14118301	5
Caelo	Anionic hydrophilic cream DAB	14118401	11
Caelo	Anionic hydrophilic cream DAB	14118501	9
Caelo	Anionic hydrophilic cream DAB	14138801	5
Caelo	Anionic hydrophilic cream DAB	14139001	2
Caelo	Anionic hydrophilic cream DAB	14138902	5
Caelo	Anionic hydrophilic cream DAB	14138901	5
Caelo	Anionic hydrophilic cream DAB	14144301	10
Caelo	Anionic hydrophilic cream DAB	14144501	5
Caelo	Anionic hydrophilic cream DAB	14146201	5
Caelo	Anionic hydrophilic cream DAB	14144701	3
Caelo	Anionic hydrophilic cream DAB	14146301	3
Caelo	Anionic hydrophilic cream DAB	14146401	3
Phönix Gotha	Hydrophilic cream, non-ionic...	141521	2
Caelo	Anionic hydrophilic cream DAB	14185401	5
Caelo	Anionic hydrophilic cream DAB	14185501	5
Caelo	Anionic hydrophilic cream DAB	14185701	1
Caelo	Anionic hydrophilic cream DAB	14185801	6
Caelo	Anionic hydrophilic cream DAB	14186601	8
Fagron	Anionic hydrophilic cream DAB	14186601	1
Phoenix 29.08.2014 5,78	Anionic hydrophilic cream DAB	14186601	1
Caelo	Anionic hydrophilic cream DAB	14186702	7
Phönix Gotha	Hydrophilic cream, non-ionic...	142311	2
Sana 08.10.14 EK 5,76	Anionic hydrophilic cream DAB	14232901	1
Caelo	Anionic hydrophilic cream DAB	14233001	4
Caelo	Anionic hydrophilic cream DAB	14233101	12
Caelo	Anionic hydrophilic cream DAB	14233201	7
Caelo	Anionic hydrophilic cream DAB	14233202	9
Caelo	Anionic hydrophilic cream DAB	14233401	10
Caelo	Anionic hydrophilic cream DAB	14233601	5
Caelo	Anionic hydrophilic cream DAB	14263901	7
Caelo	Anionic hydrophilic cream DAB	14264001	8
Caelo	Anionic hydrophilic cream DAB	14264101	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	14264201	8
Caelo	Anionic hydrophilic cream DAB	14264301	4
Anzag 27.03.2015	Anionic hydrophilic cream DAB	14264301	1
Caelo	Anionic hydrophilic cream DAB	14267401	1
Caelo	Anionic hydrophilic cream DAB	14285501	4
Caelo	Anionic hydrophilic cream DAB	14285601	7
Caelo	Anionic hydrophilic cream DAB	14286101	5
Caelo	Anionic hydrophilic cream DAB	14286001	3
Caelo	Anionic hydrophilic cream DAB	14285902	2
Caelo	Anionic hydrophilic cream DAB	14286201	8
Caelo	Anionic hydrophilic cream DAB	14294201	15
Caelo	Anionic hydrophilic cream DAB	14294501	2
Sana 11.02.15 EK 5,75	Anionic hydrophilic cream DAB	14294902	1
Caelo	Anionic hydrophilic cream DAB	14294902	2
Caelo	Anionic hydrophilic cream DAB	14295001	7
Phönix Gotha	Hydrophilic cream, non-ionic...	143011	3
Caelo	Hydrophilic cream, non-ionic...	14304902	7
Henry Lamotte Oils GmbH	Anionic hydrophilic cream DAB	14341901	1
Caelo	Anionic hydrophilic cream DAB	14341501	5
Caelo	Anionic hydrophilic cream DAB	14342201	7
Caelo	Hydrophilic cream, non-ionic...	14336001	2
Caelo	Anionic hydrophilic cream DAB	14342301	2
Caelo	Anionic hydrophilic cream DAB	14342501	6
Caelo	Anionic hydrophilic cream DAB	14342601	4
Caelo	Anionic hydrophilic cream DAB	14342701	2
Caelo	Anionic hydrophilic cream DAB	14342801	11
Caelo	Anionic hydrophilic cream DAB	14366301	5
Caelo	Anionic hydrophilic cream DAB	14366801	22
Caelo	Anionic hydrophilic cream DAB	14367001	13
Caelo	Anionic hydrophilic cream DAB	14367401	11
Phönix Gotha	Hydrophilic cream, non-ionic...	144811	5
Caelo	Anionic hydrophilic cream DAB	1485401	1
Fagron	Hydrophilic cream, non-ionic...	14D28-M01	1
Caelo	Anionic hydrophilic cream DAB	15001101	6
Caelo	Anionic hydrophilic cream DAB	15001301	10
Caelo	Anionic hydrophilic cream DAB	15001401	8
Caelo	Anionic hydrophilic cream DAB	1500141	3
Caelo	Anionic hydrophilic cream DAB	15026901	4
Caelo	Anionic hydrophilic cream DAB	15027001	27
Caelo	Anionic hydrophilic cream DAB	15027002	2
Caelo	Anionic hydrophilic cream DAB	15027101	5
Caelo	Anionic hydrophilic cream DAB	15027301	3
Caelo	Anionic hydrophilic cream DAB	15027201	4
Caelo	Anionic hydrophilic cream DAB	15034201	9
Caelo	Anionic hydrophilic cream DAB	15034301	7
Caelo	Anionic hydrophilic cream DAB	15034601	7
Caelo	Anionic hydrophilic cream DAB	15034401	5
Caelo	Hydrophilic cream, non-ionic...	15037703	8
Caelo	Hydrophilic cream, non-ionic...	15039001	1
Phönix Gotha	Hydrophilic cream, non-ionic...	150751	1
Caelo	Hydrophilic cream, non-ionic...	15039001	7
Caelo	Hydrophilic cream, non-ionic...	15155904	4
Caelo	Hydrophilic cream, non-ionic...	15155904	1
Phönix Gotha	Hydrophilic cream, non-ionic...	151711	2
Leyh Pharma / Phönix	Hydrophilic cream, non-ionic...	151711	1
Caelo	Anionic hydrophilic cream DAB	15199801	2
Caelo	Anionic hydrophilic cream DAB	15200101	8
Caelo	Anionic hydrophilic cream DAB	15200301	6
Caelo	Anionic hydrophilic cream DAB	15200201	5
Caelo	Anionic hydrophilic cream DAB	15200501	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic cream, non-ionic...	15210303	6
Caelo	Anionic hydrophilic cream DAB	15258901	4
Caelo	Anionic hydrophilic cream DAB	15259001	8
Caelo	Anionic hydrophilic cream DAB	15259101	6
Caelo	Anionic hydrophilic cream DAB	15266901	3
Caelo	Anionic hydrophilic cream DAB	15267101	3
Caelo	Anionic hydrophilic cream DAB	15267201	7
Caelo	Anionic hydrophilic cream DAB	15267301	6
Caelo	Anionic hydrophilic cream DAB	15267501	1
Caelo	Anionic hydrophilic cream DAB	15267601	3
phoenix/noveda	Anionic hydrophilic cream DAB	15267501	1
Caelo	Anionic hydrophilic cream DAB	15267801	8
Phönix Gotha	Hydrophilic cream, non-ionic...	153041	4
Caelo	Anionic hydrophilic cream DAB	15315501	3
Caelo	Anionic hydrophilic cream DAB	15315701	7
Caelo	Anionic hydrophilic cream DAB	15315801	2
Caelo	Anionic hydrophilic cream DAB	15315901	3
Caesar & Loretz/AHD	Anionic hydrophilic cream DAB	15319101	1
Caelo	Anionic hydrophilic cream DAB	15319101	3
Caelo	Anionic hydrophilic cream DAB	15319201	5
Caelo	Anionic hydrophilic cream DAB	15319301	6
Caelo	Anionic hydrophilic cream DAB	15319501	7
Caelo	Hydrophilic cream, non-ionic...	15325702	3
Caelo	Hydrophilic cream, non-ionic...	15325801	13
Caelo	Anionic hydrophilic cream DAB	15352901	9
Caelo	Anionic hydrophilic cream DAB	15353201	3
Caelo	Anionic hydrophilic cream DAB	15357301	3
Caelo	Anionic hydrophilic cream DAB	15357401	5
Caesar&Loretz GmbH	Anionic hydrophilic cream DAB	15357401	1
Phönix Gotha	Hydrophilic cream, non-ionic...	153621	1
Cael Sanacorp 05.02.2016	Anionic hydrophilic cream DAB	15357401	2
LEYH PHARMA / Phönix Le...	Hydrophilic cream, non-ionic...	153621	1
Caelo	Anionic hydrophilic cream DAB	15396801	3
Caelo	Anionic hydrophilic cream DAB	15397301	9
Caelo	Anionic hydrophilic cream DAB	15397201	5
Caelo	Anionic hydrophilic cream DAB	15400801	1
Caelo	Anionic hydrophilic cream DAB	15400401	3
Caelo	Anionic hydrophilic cream DAB	154009	1
Caelo	Anionic hydrophilic cream DAB	15400901	4
Caelo	Anionic hydrophilic cream DAB	15401001	6
Caelo	Anionic hydrophilic cream DAB	15401101	7
Caelo	Anionic hydrophilic cream DAB	15401201	1
Euro OTC	Anionic hydrophilic cream DAB	15401302	1
Caelo	Anionic hydrophilic cream DAB	15401401	13
Cealo/Noweda	Anionic hydrophilic cream DAB	15401401	2
Caelo	Anionic hydrophilic cream DAB	15401701	5
Caelo	Anionic hydrophilic cream DAB	15428102	4
Caelo	Anionic hydrophilic cream DAB	1541702	1
Caelo	Anionic hydrophilic cream DAB	15428201	8
Caelo	Anionic hydrophilic cream DAB	15428401	3
Noweda	Anionic hydrophilic cream DAB	15428401	1
Caelo	Anionic hydrophilic cream DAB	15431702	11
Caelo	Anionic hydrophilic cream DAB	15432101	6
Caelo	Anionic hydrophilic cream DAB	15432001	6
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	155221	1
Leyh Pharma / Alliance ...	Hydrophilic cream, non-ionic...	155221	1
LEYH PHARMA / Phönix Le...	Hydrophilic cream, non-ionic...	155221	1
apomix/Holdermann	Hydrophilic cream, non-ionic...	15B290615	1
Caelo	Hydrophilic cream, non-ionic...	16017501	6
Leyh Pharma/Kehr	Hydrophilic cream, non-ionic...	16017501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	16046301	3
Caelo	Anionic hydrophilic cream DAB	16046801	2
Caelo	Anionic hydrophilic cream DAB	16047101	1
Caelo	Hydrophilic cream, non-ionic...	16054101	4
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	160651	1
Leyh Pharma/Anzag	Hydrophilic cream, non-ionic...	160651	1
Caelo	Anionic hydrophilic cream DAB	16103201	1
Caelo	Hydrophilic cream, non-ionic...	16114102	4
Caelo	Hydrophilic cream, non-ionic...	16114201	6
Anzag	Anionic hydrophilic cream DAB	16158701	3
AHD	Anionic hydrophilic cream DAB	16180101	1
Caelo	Hydrophilic cream, non-ionic...	16231701	1
Caelo	Hydrophilic cream, non-ionic...	16222001	6
Leyh Pharma / Phönix	Hydrophilic cream, non-ionic...	162241	1
Caelo	Anionic hydrophilic cream DAB	1624701	1
Caelo	Hydrophilic cream, non-ionic...	16257902	5
Gehe	Anionic hydrophilic cream DAB	16250201	1
Gehe	Anionic hydrophilic cream DAB	16249902	1
Sanacorp	Anionic hydrophilic cream DAB	16258901	1
Ebert-Jacobi	Anionic hydrophilic cream DAB	16259201	1
Caelo	Anionic hydrophilic cream DAB	16259201	1
AHD	Anionic hydrophilic cream DAB	16259201	1
Caelo	Hydrophilic cream, non-ionic...	16321701	6
Caelo	Anionic hydrophilic cream DAB	16322601	1
Leyh Pharma / Phönix	Hydrophilic cream, non-ionic...	164441	1
Bombastus	Hydrophilic cream, non-ionic...	17000743	3
Tsiripidou	Hydrophilic cream, non-ionic...	17000748	1
Caelo	Hydrophilic cream, non-ionic...	17009601	1
Caelo	Hydrophilic cream, non-ionic...	1700971	1
Bombastus	Anionic hydrophilic cream DAB	17001308	1
Caelo	Anionic hydrophilic cream DAB	17030001	2
Phoenix	Anionic hydrophilic cream DAB	17054201	1
Caelo	Anionic hydrophilic cream DAB	17088501	2
Noweda	Anionic hydrophilic cream DAB	17088901	1
Caelo	Hydrophilic cream, non-ionic...	17113403	2
Caelo	Anionic hydrophilic cream DAB	17112014B	1
Caelo	Anionic hydrophilic cream DAB	17142002	1
Caelo	Anionic hydrophilic cream DAB	17142701	1
Caelo	Anionic hydrophilic cream DAB	17143201	1
Caelo	Hydrophilic cream, non-ionic...	17168502	1
AHD	Anionic hydrophilic cream DAB	17200701	2
Caelo	Anionic hydrophilic cream DAB	17201001	1
Noweda	Anionic hydrophilic cream DAB	17201501	1
Noweda	Anionic hydrophilic cream DAB	17202101	1
Leyh Pharma/Sanacorp	Hydrophilic cream, non-ionic...	172051	1
Caelo	Anionic hydrophilic cream DAB	17248601	1
GEHE	Anionic hydrophilic cream DAB	17270001	1
Noweda	Anionic hydrophilic cream DAB	17270401	1
Caelo	Anionic hydrophilic cream DAB	17299601	1
Caelo	Hydrophilic cream, non-ionic...	17313903	1
apomix	Hydrophilic cream, non-ionic...	17A100315	1
Apomix / Sanacorp	Hydrophilic cream, non-ionic...	17B040815	3
Apomix /Noweda	Hydrophilic cream, non-ionic...	17b040815	1
Gehe	Hydrophilic cream, non-ionic...	17B120814	1
GEHE	Anionic hydrophilic cream DAB	18066901	1
Apomix	Hydrophilic cream, non-ionic...	18B270916	2
Caelo	Anionic hydrophilic cream DAB	1814E2540	1
Bombastus	Hydrophilic cream, non-ionic...	201015	3
PKH GmbH	Hydrophilic cream, non-ionic...	21A110515	1
Pharmachem	Hydrophilic cream, non-ionic...	220713	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic cream, non-ionic...	22A050917	1
apomix	Hydrophilic cream, non-ionic...	25B081215	1
Apomix / Noweda	Hydrophilic cream, non-ionic...	23B101115	1
Allianc. Health. /Apomix	Hydrophilic cream, non-ionic...	25A300516	1
GEHE/Apomix	Hydrophilic cream, non-ionic...	27A220715	2
Bombastus	Anionic hydrophilic cream DAB	277463	1
Bombastus	Anionic hydrophilic cream DAB	280291	1
Bombastus	Anionic hydrophilic cream DAB	282003	1
Caelo	Anionic hydrophilic cream DAB	281999	1
Bombastus	Anionic hydrophilic cream DAB	282006	5
Bombastus	Anionic hydrophilic cream DAB	282110	5
Bombastus	Anionic hydrophilic cream DAB	282109	1
Bombastus	Anionic hydrophilic cream DAB	282114	4
Bombastus	Anionic hydrophilic cream DAB	283411	4
Raben-Apotheke Mildena	Anionic hydrophilic cream DAB	283410	1
Bombastus	Anionic hydrophilic cream DAB	283413	2
Bombastus	Anionic hydrophilic cream DAB	283474	1
Bombastus	Anionic hydrophilic cream DAB	283477	2
Bomba	Anionic hydrophilic cream DAB	283478	1
Bombastus	Anionic hydrophilic cream DAB	283482	2
Bombastus	Hydrophilic cream, non-ionic...	283638	2
Bombastus	Hydrophilic cream, non-ionic...	283695	6
Bombastus	Hydrophilic cream, non-ionic...	283639	2
Bombastus	Hydrophilic cream, non-ionic...	284899	4
Bombastus	Anionic hydrophilic cream DAB	285592	1
Bombastus	Anionic hydrophilic cream DAB	285591	1
Bombastus	Anionic hydrophilic cream DAB	285595	1
Bombastus	Anionic hydrophilic cream DAB	285596	2
Bombastus	Anionic hydrophilic cream DAB	285594	2
Bombastus	Anionic hydrophilic cream DAB	285597	2
Bombastus	Hydrophilic cream, non-ionic...	286326	7
Bombastus	Hydrophilic cream, non-ionic...	286330	5
Bombastus	Anionic hydrophilic cream DAB	286498	2
Bombastus	Anionic hydrophilic cream DAB	286495	2
Bombastus	Anionic hydrophilic cream DAB	286499	3
Bombastus	Anionic hydrophilic cream DAB	286502	3
Bombastus	Anionic hydrophilic cream DAB	286505	2
Bombastus	Anionic hydrophilic cream DAB	287432	2
Bombastus	Anionic hydrophilic cream DAB	287433	1
Bombastus	Anionic hydrophilic cream DAB	287557	6
Bombastus	Anionic hydrophilic cream DAB	287436	2
Bombastus	Anionic hydrophilic cream DAB	287559	2
Bombastus	Anionic hydrophilic cream DAB	287562	4
Caelo	Anionic hydrophilic cream DAB	287562	1
Bombastus	Anionic hydrophilic cream DAB	287563	1
Bombastus	Anionic hydrophilic cream DAB	287564	2
Bombastus	Anionic hydrophilic cream DAB	287565	1
Bombastus	Hydrophilic cream, non-ionic...	287715	3
Bombastus	Hydrophilic cream, non-ionic...	287763	3
Bombastus	Hydrophilic cream, non-ionic...	287761	4
Bombastus	Anionic hydrophilic cream DAB	288570	1
Bombastus	Anionic hydrophilic cream DAB	288573	1
Phönix	Anionic hydrophilic cream DAB	288572	1
Bombastus	Anionic hydrophilic cream DAB	288574	3
Bombastus	Anionic hydrophilic cream DAB	288652	1
Sanacorp	Anionic hydrophilic cream DAB	288573	1
Bombastus	Anionic hydrophilic cream DAB	288653	1
Bombastus	Anionic hydrophilic cream DAB	288655	2
Bombastus	Anionic hydrophilic cream DAB	288656	2
Bombastus	Anionic hydrophilic cream DAB	288658	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Anionic hydrophilic cream DAB	288658	1
Bombastus	Anionic hydrophilic cream DAB	288657	1
Bombastus	Anionic hydrophilic cream DAB	288721	2
Bombastus	Anionic hydrophilic cream DAB	288659	3
Bombastus	Hydrophilic cream, non-ionic...	288795	4
Caelo	Hydrophilic cream, non-ionic...	28882	1
Caelo	Hydrophilic cream, non-ionic...	289942	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	28A220616	1
Bombastus	Hydrophilic cream, non-ionic...	289944	1
Bombastus	Anionic hydrophilic cream DAB	290078	1
Bombastus	Anionic hydrophilic cream DAB	290079	1
Bombastus	Anionic hydrophilic cream DAB	290080	2
Bombastus	Anionic hydrophilic cream DAB	290144	2
Bombastus	Anionic hydrophilic cream DAB	290145	2
Bombastus	Anionic hydrophilic cream DAB	290146	4
Rabenapotheke	Anionic hydrophilic cream DAB	290146	1
Noweda	Anionic hydrophilic cream DAB	290147	2
Bombastus	Anionic hydrophilic cream DAB	290220	2
Bombastus	Anionic hydrophilic cream DAB	290150	1
Bombastus	Hydrophilic cream, non-ionic...	291384	3
Bombastus	Hydrophilic cream, non-ionic...	291385	3
Bombastus	Hydrophilic cream, non-ionic...	291386	3
Bombastus	Anionic hydrophilic cream DAB	291650	1
Bombastus	Anionic hydrophilic cream DAB	291656	4
Bombastus	Anionic hydrophilic cream DAB	291657	1
Bombastus	Anionic hydrophilic cream DAB	291658	4
Bobmastus	Anionic hydrophilic cream DAB	291657	1
Bombastus	Anionic hydrophilic cream DAB	291659	2
Bombastus	Anionic hydrophilic cream DAB	291731	1
Bombastus	Anionic hydrophilic cream DAB	291729	1
Bombastus	Anionic hydrophilic cream DAB	291732	2
Bombastus	Anionic hydrophilic cream DAB	291822	2
Bombastus	Anionic hydrophilic cream DAB	291823	6
Bombastus	Anionic hydrophilic cream DAB	291824	6
Bombastus	Anionic hydrophilic cream DAB	291825	3
Bombastus	Hydrophilic cream, non-ionic...	292452	2
Bombastus	Hydrophilic cream, non-ionic...	292739	2
Bombastus	Hydrophilic cream, non-ionic...	292736	1
Bombastus	Hydrophilic cream, non-ionic...	292804	6
Bombastus	Hydrophilic cream, non-ionic...	292807	6
Bombastus	Hydrophilic cream, non-ionic...	292811	5
Bombastus	Hydrophilic cream, non-ionic...	292812	3
Bombastus	Hydrophilic cream, non-ionic...	292814	11
Bombastus	Hydrophilic cream, non-ionic...	292970	1
Bombastus	Anionic hydrophilic cream DAB	293357	1
Bombastus	Anionic hydrophilic cream DAB	293359	4
Caelo	Anionic hydrophilic cream DAB	293437	1
Phoenix Goettin	Anionic hydrophilic cream DAB	293439	1
Bombastus	Anionic hydrophilic cream DAB	293440	5
Bombastus	Anionic hydrophilic cream DAB	293437	2
Bombastus	Anionic hydrophilic cream DAB	293441	6
Caelo	Anionic hydrophilic cream DAB	293440	1
Bombastus	Anionic hydrophilic cream DAB	293447	3
Bombastus	Anionic hydrophilic cream DAB	293445	1
Bombastus	Hydrophilic cream, non-ionic...	294310	2
Bombastus	Anionic hydrophilic cream DAB	293448	3
Bombastus	Hydrophilic cream, non-ionic...	294406	7
Bombastus	Hydrophilic cream, non-ionic...	294407	5
Bombastus	Anionic hydrophilic cream DAB	294545	4
Bombastus	Anionic hydrophilic cream DAB	294546	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	294547	2
Bombastus	Anionic hydrophilic cream DAB	294641	1
Bombastus	Anionic hydrophilic cream DAB	294549	2
Bombastus	Anionic hydrophilic cream DAB	294643	5
Bombastus	Anionic hydrophilic cream DAB	294642	1
Bombastus	Anionic hydrophilic cream DAB	294644	5
Caelo	Anionic hydrophilic cream DAB	294645	1
Bombastus	Anionic hydrophilic cream DAB	294647	3
Bombastus	Anionic hydrophilic cream DAB	294646	2
Bombastus	Anionic hydrophilic cream DAB	2946450	1
Bombastus	Anionic hydrophilic cream DAB	294648	1
Caelo	Anionic hydrophilic cream DAB	294649	2
Bombastus	Anionic hydrophilic cream DAB	294649	5
Bombastus	Anionic hydrophilic cream DAB	294650	7
Bombastus	Anionic hydrophilic cream DAB	295396	1
Bombastus	Anionic hydrophilic cream DAB	295397	2
Bombastus	Anionic hydrophilic cream DAB	295400	6
Caelo	Anionic hydrophilic cream DAB	295400	1
Bombastus	Anionic hydrophilic cream DAB	295403	1
Bombastus	Anionic hydrophilic cream DAB	295511	5
Bombastus	Anionic hydrophilic cream DAB	295512	7
Bombastus	Anionic hydrophilic cream DAB	295514	2
Bombastus	Anionic hydrophilic cream DAB	295515	1
Bombastus	Anionic hydrophilic cream DAB	2955126	1
Bombastus	Anionic hydrophilic cream DAB	295516	5
Bombastus	Anionic hydrophilic cream DAB	295517	8
Bombastus	Anionic hydrophilic cream DAB	295518	7
Caesar & Loretz/Noweda	Hydrophilic cream, non-ionic...	295648	1
Bombastus	Hydrophilic cream, non-ionic...	295654	4
Bombastus	Hydrophilic cream, non-ionic...	295652	2
Bombastus	Hydrophilic cream, non-ionic...	295656	2
Bombastus	Anionic hydrophilic cream DAB	296700	4
Bombastus	Hydrophilic cream, non-ionic...	295657	8
Bombastus	Anionic hydrophilic cream DAB	296699	3
Bombastus	Anionic hydrophilic cream DAB	296702	6
Bombastus	Anionic hydrophilic cream DAB	296703	8
Bombastus	Anionic hydrophilic cream DAB	296705	3
Bombastus	Anionic hydrophilic cream DAB	296707	3
Bombastus	Anionic hydrophilic cream DAB	296706	1
Bombastus	Anionic hydrophilic cream DAB	296708	3
Bombastus	Anionic hydrophilic cream DAB	296782	1
Bombastus	Anionic hydrophilic cream DAB	296783	4
Bombastus	Anionic hydrophilic cream DAB	296784	2
Bombastus	Anionic hydrophilic cream DAB	296786	2
Bombastus	Anionic hydrophilic cream DAB	296787	6
Bombastus	Hydrophilic cream, non-ionic...	297305	6
Bombastus	Hydrophilic cream, non-ionic...	297298	1
Bombastus	Hydrophilic cream, non-ionic...	297302	1
Bombastus	Hydrophilic cream, non-ionic...	297304	3
Bombastus	Hydrophilic cream, non-ionic...	297309	8
Bombastus	Hydrophilic cream, non-ionic...	297310	2
Bombastus	Hydrophilic cream, non-ionic...	297370	1
Bombastus	Hydrophilic cream, non-ionic...	297372	6
Fagron	Hydrophilic cream, non-ionic...	297372	1
Bombastus	Hydrophilic cream, non-ionic...	297373	4
Bombastus	Hydrophilic cream, non-ionic...	297374	2
Bombastus	Anionic hydrophilic cream DAB	297662	5
Bombastus	Anionic hydrophilic cream DAB	297663	1
Bombastus	Anionic hydrophilic cream DAB	297668	6
Bombastus	Anionic hydrophilic cream DAB	297671	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	297669	5
Bombastus	Anionic hydrophilic cream DAB	297672	2
Bombastus	Anionic hydrophilic cream DAB	297675	5
Bombastus	Anionic hydrophilic cream DAB	297676	2
Bombastus	Anionic hydrophilic cream DAB	297679	6
Bombastus	Anionic hydrophilic cream DAB	297677	1
Bombastus	Anionic hydrophilic cream DAB	297680	4
Bombastus	Hydrophilic cream, non-ionic...	297911	1
Bombastus	Hydrophilic cream, non-ionic...	298978	1
Bombastus	Hydrophilic cream, non-ionic...	298981	6
Bombastus	Hydrophilic cream, non-ionic...	3x297305	1
Euro OTC	Anionic hydrophilic cream DAB	300352	1
Bombastus	Hydrophilic cream, non-ionic...	300758	1
Bombastus	Hydrophilic cream, non-ionic...	300756	3
Bombastus	Hydrophilic cream, non-ionic...	300761	3
Bombastus	Hydrophilic cream, non-ionic...	300763	5
Bombastus	Hydrophilic cream, non-ionic...	300764	4
Bombastus	Hydrophilic cream, non-ionic...	300764x2	1
Bombastus	Hydrophilic cream, non-ionic...	300766	9
Bombastus	Hydrophilic cream, non-ionic...	301622	5
Bombastus	Hydrophilic cream, non-ionic...	301623	5
Bombastus	Hydrophilic cream, non-ionic...	301624	2
Bombastus	Hydrophilic cream, non-ionic...	301627	10
Bombastus	Anionic hydrophilic cream DAB	302452	3
Bombastus	Anionic hydrophilic cream DAB	302460	1
Bombastus	Hydrophilic cream, non-ionic...	302690	5
Bombastus	Hydrophilic cream, non-ionic...	302684	3
Bombastus	Hydrophilic cream, non-ionic...	302686	2
Bombastus	Hydrophilic cream, non-ionic...	302688	6
Bombastus	Hydrophilic cream, non-ionic...	303735	1
Bombastus	Hydrophilic cream, non-ionic...	304771	1
Bombastus	Hydrophilic cream, non-ionic...	304773	1
Bombastus	Hydrophilic cream, non-ionic...	304774	7
Bombastus	Hydrophilic cream, non-ionic...	304775	4
Bombastus	Hydrophilic cream, non-ionic...	304776	3
Bombastus	Anionic hydrophilic cream DAB	305370	5
Bombastus	Hydrophilic cream, non-ionic...	304777	3
Bombastus	Anionic hydrophilic cream DAB	305377	1
Bombastus	Anionic hydrophilic cream DAB	305638	1
Bombastus	Hydrophilic cream, non-ionic...	306126	5
Bombastus	Hydrophilic cream, non-ionic...	306123	1
Bombastus	Hydrophilic cream, non-ionic...	306127	2
Bombastus	Hydrophilic cream, non-ionic...	306130	2
Bombastus/San	Hydrophilic cream, non-ionic...	306131	1
Bombastus	Hydrophilic cream, non-ionic...	306131	10
Bombastus	Anionic hydrophilic cream DAB	307024	1
Bombastus	Hydrophilic cream, non-ionic...	307531	7
Bombastus	Hydrophilic cream, non-ionic...	307526	2
Bombastus	Anionic hydrophilic cream DAB	308392	1
Bombastus	Hydrophilic cream, non-ionic...	307530	3
Bombastus	Anionic hydrophilic cream DAB	308396	3
Bombastus	Anionic hydrophilic cream DAB	308398	2
Bombastus	Anionic hydrophilic cream DAB	308404	1
Bombastus	Hydrophilic cream, non-ionic...	309904	5
Bombastus	Hydrophilic cream, non-ionic...	309908	6
Bombastus	Anionic hydrophilic cream DAB	310266	3
Bombastus	Hydrophilic cream, non-ionic...	310515	2
Bombastus	Hydrophilic cream, non-ionic...	310518	4
Bombastus	Hydrophilic cream, non-ionic...	310520	6
Bombastus	Anionic hydrophilic cream DAB	311123	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix / Anzag	Hydrophilic cream, non-ionic...	32A250414	1
APOMIX/AHD	Hydrophilic cream, non-ionic...	35A130514	1
Caelo	Hydrophilic cream, non-ionic...	36A140514	1
apomix / Anzag	Hydrophilic cream, non-ionic...	36A140514	1
Apomix/Kehr	Hydrophilic cream, non-ionic...	36A261016	1
Fagron	Hydrophilic cream, non-ionic...	39131216	1
Caelo	Hydrophilic cream, non-ionic...	47A071014	1
apomix / Anzag	Hydrophilic cream, non-ionic...	51A261114	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	47A2801015	1
Allianc. Health. / Apomix	Hydrophilic cream, non-ionic...	48A181115	1
Bombastus	Anionic hydrophilic cream DAB	60080602	2
Bombastus	Hydrophilic cream, non-ionic...	60089512	1
Caelo	Anionic hydrophilic cream DAB	9/AKD13	1
Sanacorp	Anionic hydrophilic cream DAB	AKD11	1
Bombastus	Hydrophilic cream, non-ionic...	99320	1
Caelo	Anionic hydrophilic cream DAB	AKD10	1
Bombastus	Hydrophilic cream, non-ionic...	268347	1
Bombastus	Hydrophilic cream, non-ionic...	269714	2
Gehe	Hydrophilic cream, non-ionic...	271518	1
Bombastus	Hydrophilic cream, non-ionic...	271522	2
Bombastus	Hydrophilic cream, non-ionic...	272492	2
Bombastus	Hydrophilic cream, non-ionic...	273130	4
Bombastus	Hydrophilic cream, non-ionic...	273134	2
Bombastus	Hydrophilic cream, non-ionic...	274219	2
Bombastus	Hydrophilic cream, non-ionic...	274103	2
Bombastus	Hydrophilic cream, non-ionic...	274216	3
GEHE/apomix/22.05.2012	Hydrophilic cream, non-ionic...	3A130212	1
Bombastus	Hydrophilic cream, non-ionic...	274223	4
Bombastus	Hydrophilic cream, non-ionic...	274211	5
Phönix Gotha	Hydrophilic cream, non-ionic...	4B070213	1
Pharmachem GmbH	Hydrophilic cream, non-ionic...	50912	1
Apomix	Hydrophilic cream, non-ionic...	8B150413	1
Noweda	Hydrophilic cream, non-ionic...	11A260312	1
Gehe	Hydrophilic cream, non-ionic...	121431	1
Leyh Pharma	Hydrophilic cream, non-ionic...	122841	1
Phönix Gotha	Hydrophilic cream, non-ionic...	122841	4
Phönix	Hydrophilic cream, non-ionic...	124511	1
Phönix Gotha	Hydrophilic cream, non-ionic...	124511	5
Leyh Pharma	Hydrophilic cream, non-ionic...	130651	2
Phönix Gotha	Hydrophilic cream, non-ionic...	133121	3
Phönix Gotha	Hydrophilic cream, non-ionic...	131451	2
Phönix Gotha	Hydrophilic cream, non-ionic...	134841	2
Fagron	Hydrophilic cream, non-ionic...	13K12-M04	1
Fagron	Hydrophilic cream, non-ionic...	13A250313	1
Fagron	Hydrophilic cream, non-ionic...	13D05-M08	2
Fagron	Hydrophilic cream, non-ionic...	13K12-M05	1
GEHE/apomix/13.07.2012	Hydrophilic cream, non-ionic...	13A230412	1
Gehe	Hydrophilic cream, non-ionic...	14B090713	1
Apomix	Hydrophilic cream, non-ionic...	16A300512	2
Apomix	Hydrophilic cream, non-ionic...	18A070612	1
Pharmachem GmbH	Hydrophilic cream, non-ionic...	220713	2
Caelo	Hydrophilic cream, non-ionic...	274223	1
Bombastus	Hydrophilic cream, non-ionic...	276252	7
Caelo	Hydrophilic cream, non-ionic...	276596	1
Bombastus	Hydrophilic cream, non-ionic...	277873	4
Bombastus	Hydrophilic cream, non-ionic...	277637	1
Bombastus	Hydrophilic cream, non-ionic...	277879	1
Bombastus	Hydrophilic cream, non-ionic...	277878	11
Bombastus	Hydrophilic cream, non-ionic...	278027	3
Bombastus	Hydrophilic cream, non-ionic...	279661	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrophilic cream, non-ionic...	279655	2
Apomix	Hydrophilic cream, non-ionic...	27B221012	1
Caelo	Hydrophilic cream, non-ionic...	280988	1
Bombastus	Hydrophilic cream, non-ionic...	282216	8
Gehe	Hydrophilic cream, non-ionic...	282214	1
gehe	Hydrophilic cream, non-ionic...	282211	1
Bombastus	Hydrophilic cream, non-ionic...	283636	6
Bombastus	Hydrophilic cream, non-ionic...	283692	2
Bombastus	Hydrophilic cream, non-ionic...	2886797	1
Bombastus	Hydrophilic cream, non-ionic...	283637	2
Apomix / Phönix	Hydrophilic cream, non-ionic...	32A090713	1
Apomix	Hydrophilic cream, non-ionic...	39A131112	1
apomix/Anzag	Hydrophilic cream, non-ionic...	40A281013	2
Apomix	Hydrophilic cream, non-ionic...	36B270513	1
Fagron	Hydrophilic cream, non-ionic...	4A220113	1
Apomix	Hydrophilic cream, non-ionic...	53A191113	1
Caelo	Anionic hydrophilic cream DAB	2569516	1
Bombastus	Anionic hydrophilic cream DAB	263981	5
Bombastus	Anionic hydrophilic cream DAB	262774	1
Bombastus	Anionic hydrophilic cream SR...	264643	1
Bombastus	Anionic hydrophilic cream DAB	264833	1
Klenk	Anionic hydrophilic cream DAB	264833	1
Bombastus	Anionic hydrophilic cream SR...	264644	1
Klenk	Anionic hydrophilic cream DAB	266259	1
Bombastus	Anionic hydrophilic cream DAB	267616	1
Bombastus	Anionic hydrophilic cream DAB	267436	1
Bombastus	Anionic hydrophilic cream DAB	267442	1
Gehe	Anionic hydrophilic cream SR...	267647	1
Gehe	Anionic hydrophilic cream SR...	269235	1
Bombastus	Anionic hydrophilic cream DAB	269444	1
Bombastus	Anionic hydrophilic cream DAB	269513	2
Bombastus	Anionic hydrophilic cream DAB	269518	1
Bombastus	Anionic hydrophilic cream DAB	269517	1
Bombastus	Anionic hydrophilic cream DAB	270804	1
Bombastus	Anionic hydrophilic cream DAB	270824	3
Bombastus	Anionic hydrophilic cream DAB	270947	1
DAB	Anionic hydrophilic cream DAB	270951	1
Bombastus	Anionic hydrophilic cream DAB	270950	1
Bombastus	Anionic hydrophilic cream DAB	270946	1
Bombastus	Anionic hydrophilic cream DAB	271008	2
Gehe	Anionic hydrophilic cream DAB	270949	1
Bombastus	Anionic hydrophilic cream DAB	271009	1
Gehe	Anionic hydrophilic cream DAB	271010	1
Bombastus	Anionic hydrophilic cream DAB	271011	1
Bombastus	Anionic hydrophilic cream SR...	271710	2
Bombastus	Anionic hydrophilic cream SR...	271702	1
Bombastus	Anionic hydrophilic cream DAB	273536	1
Bombastus	Anionic hydrophilic cream DAB	273407	1
Bomastus	Anionic hydrophilic cream DAB	273404	1
Bombastus	Anionic hydrophilic cream DAB	273534	1
Bombastus	Anionic hydrophilic cream DAB	273532	1
Sanacorp	Anionic hydrophilic cream SR...	273693	1
Caelo	Anionic hydrophilic cream DAB	1061302	1
Caelo	Anionic hydrophilic cream DAB	1119720	1
Caelo	Anionic hydrophilic cream DAB	204312540	1
Caelo	Anionic hydrophilic cream DAB	3061302	1
Caelo	Anionic hydrophilic cream DAB	4131I-02540	1
Caelo	Anionic hydrophilic cream DAB	6121304	1
GEHE/AMH/22.05.2012	Anionic hydrophilic cream SR...	5A200212	1
Caelo	Anionic hydrophilic cream DAB	413I-02540	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream SR. ...	1101212540	1
Caelo	Anionic hydrophilic cream DAB	11180201	1
Caelo	Anionic hydrophilic cream DAB	11120201	1
Caelo	Anionic hydrophilic cream DAB	1119387	1
Caelo	Anionic hydrophilic cream DAB	11229001	2
Noweda	Anionic hydrophilic cream DAB	11229001	2
Caelo	Anionic hydrophilic cream DAB	11229101	1
Caelo	Anionic hydrophilic cream DAB	11253401	1
Caelo	Anionic hydrophilic cream DAB	11258301	1
Caelo	Anionic hydrophilic cream DAB	11253501	1
Caelo	Anionic hydrophilic cream DAB	11258201	2
Caelo	Anionic hydrophilic cream DAB	11229201	1
Caesar & Loretz GmbH, 4. ...	Anionic hydrophilic cream DAB	11311101	3
Caesar & Loretz GmbH 40. ...	Anionic hydrophilic cream DAB	11253201	1
Caelo	Anionic hydrophilic cream DAB	11311101	2
Gehe	Anionic hydrophilic cream DAB	11311202	1
Caelo	Anionic hydrophilic cream DAB	11311501	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	11320801	1
Caelo	Anionic hydrophilic cream DAB	11320901	1
Caelo	Anionic hydrophilic cream DAB	11321101	1
Caelo	Anionic hydrophilic cream DAB	11326701	3
Caesar & Loretz GmbH D-...	Anionic hydrophilic cream DAB	11326701	1
Caelo	Anionic hydrophilic cream DAB	11331501	2
Caelo	Anionic hydrophilic cream DAB	11331601	1
Caelo	Anionic hydrophilic cream DAB	11332001	1
Caelo	Anionic hydrophilic cream DAB	11332101	1
Caelo	Anionic hydrophilic cream DAB	11332401	1
Kehr	Anionic hydrophilic cream DAB	11332401	1
Caelo	Anionic hydrophilic cream DAB	11332501	1
Caelo	Anionic hydrophilic cream DAB	113332301	1
Caelo	Anionic hydrophilic cream DAB	12010201	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	12010501	1
Caelo	Anionic hydrophilic cream DAB	12010901	1
Noweda	Anionic hydrophilic cream DAB	12010901	1
Phoenix	Anionic hydrophilic cream DAB	12011001	1
Caelo	Anionic hydrophilic cream DAB	12010801	2
Caelo	Anionic hydrophilic cream DAB	12011302	2
Caelo	Anionic hydrophilic cream DAB	12012101	1
Gehe	Anionic hydrophilic cream DAB	12026001	1
Caelo	Anionic hydrophilic cream DAB	12081301	1
Caelo	Anionic hydrophilic cream DAB	12012201	1
Caelo	Anionic hydrophilic cream DAB	12083702	1
Caelo	Anionic hydrophilic cream DAB	12098801	1
Caelo	Anionic hydrophilic cream DAB	12083802	1
Noweda	Anionic hydrophilic cream DAB	12098801	1
Kehr	Anionic hydrophilic cream DAB	12098901	1
Caelo	Anionic hydrophilic cream DAB	12099001	2
Caelo	Anionic hydrophilic cream DAB	12099101	1
Caelo	Anionic hydrophilic cream DAB	12099201	3
Caelo	Anionic hydrophilic cream DAB	12099301	1
Caelo	Anionic hydrophilic cream DAB	12099501	1
DAB 2008	Anionic hydrophilic cream DAB	12099901	1
Noweda	Anionic hydrophilic cream DAB	12099901	2
Caelo	Anionic hydrophilic cream DAB	12100201	1
Caelo	Anionic hydrophilic cream DAB	12100401	3
Noweda	Anionic hydrophilic cream DAB	12100401	1
Caelo	Anionic hydrophilic cream DAB	12100601	5
Caelo	Anionic hydrophilic cream DAB	12100501	2
Caelo	Anionic hydrophilic cream DAB	12100701	3
Noweda	Anionic hydrophilic cream DAB	12100801	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	12100901	6
Noweda	Anionic hydrophilic cream DAB	12101302	5
Caelo	Anionic hydrophilic cream DAB	12128502	4
Caelo	Anionic hydrophilic cream DAB	12101201	1
Phönix Gotha	Anionic hydrophilic cream SR. . .	121431	2
Caelo	Anionic hydrophilic cream DAB	12187701	3
Caelo	Anionic hydrophilic cream DAB	12187901	8
Caelo	Anionic hydrophilic cream DAB	12195201	8
Caesar & Loretz GmbH, 4. . .	Anionic hydrophilic cream DAB	12195201	2
Caelo	Anionic hydrophilic cream DAB	12195501	2
Caelo	Anionic hydrophilic cream DAB	12195401	4
Caelo	Anionic hydrophilic cream DAB	12204301	3
Caelo	Anionic hydrophilic cream DAB	12195601	7
Caelo	Anionic hydrophilic cream DAB	12204501	2
Großhandel	Anionic hydrophilic cream DAB	12204501	1
Caelo	Anionic hydrophilic cream DAB	12204601	1
Caelo	Anionic hydrophilic cream DAB	12218501	2
Anzag,23.10.2012	Anionic hydrophilic cream DAB	12209101	1
Caelo	Anionic hydrophilic cream DAB	12218401	2
Caelo	Anionic hydrophilic cream DAB	12218301	2
Caelo	Anionic hydrophilic cream DAB	12218601	2
Caelo	Anionic hydrophilic cream DAB	12245401	4
Caelo	Anionic hydrophilic cream DAB	12245501	1
Caelo	Anionic hydrophilic cream DAB	12245601	1
Caelo	Anionic hydrophilic cream DAB	12245701	3
Sanacorp	Anionic hydrophilic cream DAB	12245701	1
Caesar & Loretz GmbH, 4. . .	Anionic hydrophilic cream DAB	12264501	1
Caelo	Anionic hydrophilic cream DAB	12264601	1
Caelo	Anionic hydrophilic cream DAB	12264801	3
Caelo	Anionic hydrophilic cream DAB	12265301	4
Caelo	Anionic hydrophilic cream DAB	12265201	6
Caelo	Anionic hydrophilic cream DAB	12265401	3
Noweda	Anionic hydrophilic cream DAB	12264901	1
Caelo	Anionic hydrophilic cream DAB	12265501	4
Caelo	Anionic hydrophilic cream DAB	12265601	1
Caelo	Anionic hydrophilic cream DAB	12265801	2
Caelo	Anionic hydrophilic cream DAB	12265701	4
Noweda	Anionic hydrophilic cream DAB	12265701	1
Caelo	Anionic hydrophilic cream DAB	12266101	5
Caesar & Loretz GmbH, 4. . .	Anionic hydrophilic cream DAB	12313202	2
Caelo	Anionic hydrophilic cream DAB	12266201	3
Anzag 25.01.2013	Anionic hydrophilic cream DAB	12266101	1
Noweda	Anionic hydrophilic cream DAB	12313202	1
Caelo	Anionic hydrophilic cream DAB	12313202	9
Caelo	Anionic hydrophilic cream DAB	12313601	6
Caelo	Anionic hydrophilic cream DAB	12313701	8
Phoenix	Anionic hydrophilic cream DAB	12313301	1
Caelo	Anionic hydrophilic cream DAB	12313901	10
Caelo	Anionic hydrophilic cream DAB	12314001	1
Caelo	Anionic hydrophilic cream DAB	12314101	9
Anzag	Anionic hydrophilic cream DAB	12314101	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	12314201	1
GEHE	Anionic hydrophilic cream DAB	12314101	1
Caelo	Anionic hydrophilic cream DAB	12314301	5
Caelo	Anionic hydrophilic cream DAB	12314401	9
Noweda	Anionic hydrophilic cream DAB	12314601	1
Caelo	Anionic hydrophilic cream DAB	12314601	9
Caelo	Anionic hydrophilic cream DAB	12314701	8
AHD,15.04.2013	Anionic hydrophilic cream DAB	12314701	1
Hirsch	Anionic hydrophilic cream DAB	12314701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	12314801	8
Sanacorp	Anionic hydrophilic cream DAB	12314701	2
Noweda	Anionic hydrophilic cream DAB	12314801	1
Caelo	Anionic hydrophilic cream DAB	12314901	5
Caelo	Anionic hydrophilic cream DAB	12316001	1
Caelo	Anionic hydrophilic cream DAB	1303312540	1
Noweda	Anionic hydrophilic cream DAB	13002705	1
Caelo	Anionic hydrophilic cream DAB	13036001	6
Caelo	Anionic hydrophilic cream DAB	13036301	8
AHC 28.05.13 EK2,95	Anionic hydrophilic cream DAB	13036301	1
Caelo	Anionic hydrophilic cream DAB	13035901	1
Caelo	Anionic hydrophilic cream DAB	13038301	7
Caelo	Anionic hydrophilic cream DAB	13068504	1
Caelo	Anionic hydrophilic cream DAB	13070301	6
Caelo	Anionic hydrophilic cream DAB	13070401	2
Noweda	Anionic hydrophilic cream DAB	13070301	1
Caelo	Anionic hydrophilic cream DAB	13070601	7
Caelo	Anionic hydrophilic cream DAB	13070502	4
Caelo	Anionic hydrophilic cream SR...	13071202	1
Caelo	Anionic hydrophilic cream DAB	13072501	5
Caelo	Anionic hydrophilic cream DAB	13128601	7
Caelo	Anionic hydrophilic cream DAB	13128701	11
Caelo	Anionic hydrophilic cream DAB	13154001	1
Caelo	Anionic hydrophilic cream DAB	1313314401	1
Caelo	Anionic hydrophilic cream DAB	13154201	2
Caelo	Anionic hydrophilic cream DAB	13154101	4
Caelo	Anionic hydrophilic cream DAB	13154202	2
Noweda	Anionic hydrophilic cream DAB	13154301	1
Phönix	Anionic hydrophilic cream DAB	13164901	1
Caelo	Anionic hydrophilic cream DAB	13164901	4
Caelo	Anionic hydrophilic cream DAB	13165001	9
Caelo	Anionic hydrophilic cream DAB	13165601	2
Caelo	Anionic hydrophilic cream DAB	13165701	6
Caelo	Anionic hydrophilic cream DAB	13165901	7
Caelo	Anionic hydrophilic cream DAB	13165101	3
Caelo	Anionic hydrophilic cream DAB	13165801	8
Caelo	Anionic hydrophilic cream DAB	13166001	4
Caelo	Anionic hydrophilic cream DAB	13167501	12
Sanacorp	Anionic hydrophilic cream DAB	13167501	1
Caelo	Anionic hydrophilic cream DAB	13167601	2
Caelo	Anionic hydrophilic cream DAB	1316801	1
Caelo	Anionic hydrophilic cream DAB	13168101	1
Caelo	Anionic hydrophilic cream DAB	13168401	9
Caelo	Anionic hydrophilic cream DAB	13167901	1
Caelo	Anionic hydrophilic cream DAB	13168601	13
Caelo	Anionic hydrophilic cream DAB	13168801	4
Großhandel	Anionic hydrophilic cream DAB	13168801	1
Caelo	Anionic hydrophilic cream DAB	13168901	8
Caelo	Anionic hydrophilic cream DAB	13169001	10
Caelo	Anionic hydrophilic cream DAB	13169301	7
Caelo	Anionic hydrophilic cream DAB	13270701	4
Caelo	Anionic hydrophilic cream DAB	13265401	7
Caelo	Anionic hydrophilic cream DAB	13265301	8
Caelo	Anionic hydrophilic cream DAB	13271101	4
Caelo	Anionic hydrophilic cream DAB	13275601	2
Caelo	Anionic hydrophilic cream DAB	13277401	10
Caelo	Anionic hydrophilic cream DAB	13277501	4
Caelo	Anionic hydrophilic cream DAB	13280101	19
Caelo	Anionic hydrophilic cream DAB	13280201	13
Caelo	Anionic hydrophilic cream DAB	13313001	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	13312801	5
Caelo	Anionic hydrophilic cream DAB	13313201	8
Caelo	Anionic hydrophilic cream DAB	13313101	4
Caelo	Anionic hydrophilic cream DAB	13313301	5
Caelo	Anionic hydrophilic cream DAB	13313501	5
Phoenix 09.01.2014	Anionic hydrophilic cream DAB	13313501	1
Caelo	Anionic hydrophilic cream DAB	13313601	5
Caelo	Anionic hydrophilic cream DAB	13313901	10
Caelo	Anionic hydrophilic cream DAB	13313701	6
Caelo	Anionic hydrophilic cream DAB	13314001	7
Caelo	Anionic hydrophilic cream DAB	1331401	1
Caelo	Anionic hydrophilic cream DAB	13314401	10
Caelo	Anionic hydrophilic cream DAB	13314801	2
Caelo	Anionic hydrophilic cream DAB	13346601	9
Caelo	Anionic hydrophilic cream DAB	13346301	10
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	13364301	1
Caelo	Anionic hydrophilic cream DAB	13364401	5
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	14016901	2
Caelo	Anionic hydrophilic cream DAB	14062401	4
Caelo	Anionic hydrophilic cream DAB	15011301	1
Caelo	Anionic hydrophilic cream DAB	15200601	4
Caelo	Anionic hydrophilic cream SR...	1603212540	1
GEHE/apomix/05.07.2012	Anionic hydrophilic cream SR...	16A050612	1
Caelo	Anionic hydrophilic cream DAB	1913Q-02540	1
Caelo	Anionic hydrophilic cream DAB	1913R-02540	1
Pharmachem	Anionic hydrophilic cream SR...	170412	2
Caelo	Anionic hydrophilic cream DAB	2103312540	1
Caelo	Anionic hydrophilic cream DAB	2112212540	1
Caelo	Anionic hydrophilic cream DAB	2405212540	1
Bombastus	Anionic hydrophilic cream DAB	262772	1
Bombastus	Anionic hydrophilic cream DAB	266175	1
Bombastus	Anionic hydrophilic cream DAB	270704	1
Phönix	Anionic hydrophilic cream DAB	270703	1
Phoenix	Anionic hydrophilic cream DAB	270705	1
Caelo	Anionic hydrophilic cream DAB	2704212540	1
Bombastus	Anionic hydrophilic cream DAB	270706	1
Bombastus	Anionic hydrophilic cream DAB	270711	2
Caelo	Anionic hydrophilic cream DAB	27091305	1
Bombastus	Anionic hydrophilic cream DAB	274554	1
Bombastus	Anionic hydrophilic cream DAB	273401	1
Bombastus	Anionic hydrophilic cream DAB	274592	1
Bombastus	Anionic hydrophilic cream DAB	274593	2
Bombastus	Anionic hydrophilic cream DAB	274594	1
Gehe/ 28.12.12 Bombastus	Anionic hydrophilic cream DAB	274594	1
Bombastus	Anionic hydrophilic cream DAB	274595	2
Bombastus	Anionic hydrophilic cream DAB	274596	1
Bombastus	Anionic hydrophilic cream DAB	274600	2
Bombastus	Anionic hydrophilic cream DAB	274597	1
Bombastus	Anionic hydrophilic cream DAB	274695	2
Bombastus	Anionic hydrophilic cream DAB	274695/273533	1
Bombastus	Anionic hydrophilic cream DAB	274697	3
Noweda	Anionic hydrophilic cream DAB	274697	1
Bombastus	Anionic hydrophilic cream DAB	274765	2
Bombastus	Anionic hydrophilic cream DAB	274766	2
Bombastus	Anionic hydrophilic cream DAB	276320	2
Bombastus	Anionic hydrophilic cream DAB	276315	3
Bombastus	Anionic hydrophilic cream DAB	276318	1
Caelo	Anionic hydrophilic cream DAB	276321	1
Bombastus	Anionic hydrophilic cream DAB	276322	1
Bombastus	Anionic hydrophilic cream DAB	276331	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	277293	5
Bombastus	Anionic hydrophilic cream DAB	276433	1
Bombastus	Anionic hydrophilic cream DAB	277297	2
Bombastus	Anionic hydrophilic cream DAB	277295	1
Gehe 14.05.2013	Anionic hydrophilic cream DAB	277297	1
Bombastus	Anionic hydrophilic cream DAB	277298	2
Bombastus	Anionic hydrophilic cream DAB	277455	2
Bombastus	Anionic hydrophilic cream DAB	277457	6
Bombastus	Anionic hydrophilic cream DAB	277458	3
Bombastus	Anionic hydrophilic cream DAB	277461	3
Bombastus	Anionic hydrophilic cream DAB	277459	2
Caelo	Anionic hydrophilic cream DAB	277461	1
Bombastus	Anionic hydrophilic cream DAB	277465	2
Bombastus	Anionic hydrophilic cream DAB	278846	1
Caelo	Anionic hydrophilic cream DAB	278848	3
Caelo	Anionic hydrophilic cream DAB	278851	1
Bombastus	Anionic hydrophilic cream DAB	278852	3
Bombastus	Anionic hydrophilic cream DAB	279054	2
Bombastus	Anionic hydrophilic cream DAB	279055	7
Bombastus	Anionic hydrophilic cream DAB	279056	2
Bombastus	Anionic hydrophilic cream DAB	279057	1
Bombastus	Anionic hydrophilic cream DAB	279059	1
Bombastus	Anionic hydrophilic cream DAB	280289	4
Bombastus	Anionic hydrophilic cream DAB	280290	1
Bombastus	Anionic hydrophilic cream DAB	280294	1
Caelo	Anionic hydrophilic cream DAB	28031319	1
Bombastus	Anionic hydrophilic cream DAB	280418	2
Bombastus	Anionic hydrophilic cream DAB	280417	2
Bombastus	Anionic hydrophilic cream DAB	282001	1
Caelo	Anionic hydrophilic cream DAB	282005	1
Bombastus	Anionic hydrophilic cream DAB	281999	1
Gehe	Anionic hydrophilic cream DAB	282001	1
Bombastus	Anionic hydrophilic cream DAB	282004	3
Gehe	Anionic hydrophilic cream DAB	280418	2
Bombastus	Anionic hydrophilic cream DAB	282002	1
Bombastus	Anionic hydrophilic cream DAB	282005	3
Bombastus	Anionic hydrophilic cream DAB	282111	4
Bombastus	Anionic hydrophilic cream DAB	282112	1
Bombastus	Anionic hydrophilic cream DAB	282113	2
Bombastus	Anionic hydrophilic cream DAB	282115	1
Bombastus	Anionic hydrophilic cream DAB	283410	3
Bombastus	Anionic hydrophilic cream DAB	283476	1
Bombastus	Anionic hydrophilic cream DAB	283478	1
Bombastus	Anionic hydrophilic cream DAB	283481	5
Phoenix	Anionic hydrophilic cream DAB	286498	1
Caelo	Anionic hydrophilic cream SR...	2906212540	1
Caelo	Anionic hydrophilic cream DAB	30081303	1
Caelo	Anionic hydrophilic cream DAB	4313I-02540	1
Phönix Gotha	Anionic hydrophilic cream SR...	4580690	1
Caelo	Anionic hydrophilic cream SR...	3101212540	1
Caelo	Anionic hydrophilic cream DAB	654879	1
Bombastus	Anionic hydrophilic cream DAB	6861306	1
Fagron	Base cream, hydrophobic DAC	26887	1
Caelo	Base cream, hydrophobic DAC	10032015B	1
Caelo	Base cream, hydrophobic DAC	10022015F	1
Caelo	Base cream, hydrophobic DAC	1233	1
Caelo	Base cream, hydrophobic DAC	123588	1
Caelo	Base cream, hydrophobic DAC	13176002	2
Fagron	Base cream, hydrophobic DAC	13E13-T06-006	2
Fagron	Base cream, hydrophobic DAC	13T13-T06-006166	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Base cream, hydrophobic DAC	14000601	8
Hedinger/Noweda	Base cream, hydrophobic DAC	1401601	1
Caelo	Base cream, hydrophobic DAC	14028002	3
Noweda	Base cream, hydrophobic DAC	14028101	1
Caelo	Base cream, hydrophobic DAC	1406140	1
PHOENIX 11.09.2014	Base cream, hydrophobic DAC	14061501	1
caesar & Loretz GmbH	Base cream, hydrophobic DAC	14061601	2
Pharmachem GmbH	Base cream, hydrophobic DAC	140885	1
Caelo	Base cream, hydrophobic DAC	140885	1
Pharmachem	Base cream, hydrophobic DAC	140885	1
Pharmachem/ /Sanacorp	Base cream, hydrophobic DAC	140885	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	140886	1
Pharmachem	Base cream, hydrophobic DAC	140886	2
Pharmachem.Pöckneck/Phoe...	Base cream, hydrophobic DAC	140886	1
Pharmachem GmbH	Base cream, hydrophobic DAC	140886	1
Pharmachem.Pöckneck/Phoe...	Base cream, hydrophobic DAC	141318	1
PharmaChem	Base cream, hydrophobic DAC	141318	2
Pharmachem/Gehe	Base cream, hydrophobic DAC	141317	1
	Base cream, hydrophobic DAC	141318	1
Pharmachen	Base cream, hydrophobic DAC	141318	1
Nachprüfung 3083	Base cream, hydrophobic DAC	141318	1
Pharmachen GmbH/Gehe	Base cream, hydrophobic DAC	141318	1
Fischar	Base cream, hydrophobic DAC	14145301	3
Phönix	Base cream, hydrophobic DAC	14145301	1
GEHE 08.12.2014	Base cream, hydrophobic DAC	14145301	1
Caesar&Loretz GmbH/Anzag	Base cream, hydrophobic DAC	14145301	1
Caelo	Base cream, hydrophobic DAC	14145402	1
Pharmachem/ Noweda	Base cream, hydrophobic DAC	141489	1
Pharmachem/GEHE	Base cream, hydrophobic DAC	141489	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	141489	1
Pharmachen	Base cream, hydrophobic DAC	142154	3
PHARMACHEM GMBH/ALLIANC...	Base cream, hydrophobic DAC	142154	1
Pröckneck	Base cream, hydrophobic DAC	142154	2
Pharmachem/Noweda	Base cream, hydrophobic DAC	142154	2
Sanacorp/ Pharmachem	Base cream, hydrophobic DAC	142155	1
Caelo	Base cream, hydrophobic DAC	14277401	1
Bombastus	Base cream, hydrophobic DAC	143093	1
Caelo	Base cream, hydrophobic DAC	143093	1
Noweda/Pharmachem0	Base cream, hydrophobic DAC	143094	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	143093	1
Pharmachem/ Noweda	Base cream, hydrophobic DAC	143094	1
Pharmachem Pöckneck	Base cream, hydrophobic DAC	143094	2
pharma/phönix	Base cream, hydrophobic DAC	143094	1
pharmachem Pöckneck	Base cream, hydrophobic DAC	143119	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	143094	1
Pharmachem.Pöckneck/Phoe...	Base cream, hydrophobic DAC	143119	1
Caelo	Base cream, hydrophobic DAC	1435501	1
Faoron/anzag	Base cream, hydrophobic DAC	14A06-t08-009037	1
Fagron	Base cream, hydrophobic DAC	14E06-T13-012209	1
Fagron	Base cream, hydrophobic DAC	14H04-T05-014482	6
Fagron	Base cream, hydrophobic DAC	15/13-T01-026887	1
Fagron	Base cream, hydrophobic DAC	14118-T21-018124	1
Caelo	Base cream, hydrophobic DAC	14L18-T21-018124	1
Caelo	Base cream, hydrophobic DAC	15/13-T01-026887	1
Phoenix	Base cream, hydrophobic DAC	15001001	1
Caesar/Noweda	Base cream, hydrophobic DAC	15001002	1
Gehe,27.0415	Base cream, hydrophobic DAC	15001003	2
Pharmachem/Gehe	Base cream, hydrophobic DAC	150404	1
AHD/Pharmachem	Base cream, hydrophobic DAC	150404	1
Pharmachem GmbH/Noweda	Base cream, hydrophobic DAC	150405	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Pharmachem/ Gehe	Base cream, hydrophobic DAC	150404	1
Pharmachem GmbH / Anzag	Base cream, hydrophobic DAC	150405	1
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	15046401	1
Caelo	Base cream, hydrophobic DAC	1504602	1
Phoenix goettin	Base cream, hydrophobic DAC	15048701	2
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	15111402	1
Caelo	Base cream, hydrophobic DAC	15111403	3
Bombastus	Base cream, hydrophobic DAC	15111403	2
AHD /Pharmachem Pößneck	Base cream, hydrophobic DAC	151257	1
Pharmachen Pößneck	Base cream, hydrophobic DAC	151258	1
Pharamachem/Alliance He...	Base cream, hydrophobic DAC	151258	1
Pharmachem	Base cream, hydrophobic DAC	151258	2
Pharmachem GmbH	Base cream, hydrophobic DAC	151259	1
Pharmachem./ Phoenix	Base cream, hydrophobic DAC	151259	1
Caelo	Base cream, hydrophobic DAC	15132503	3
Pharmachem /Phönix	Base cream, hydrophobic DAC	152209	1
Pharmachem GmbH / Anzag	Base cream, hydrophobic DAC	152209	1
Pharmachem/ Alliance He...	Base cream, hydrophobic DAC	152209	1
AHD/Pharmachem	Base cream, hydrophobic DAC	152210	1
Pharmachem/ Gehe	Base cream, hydrophobic DAC	152210	1
Pößneck Pharma / Gehe	Base cream, hydrophobic DAC	152210	2
Caesar&Loretz GmbH /Anzag	Base cream, hydrophobic DAC	15221301	1
Caesar&Loretz GmbH	Base cream, hydrophobic DAC	15221301	1
Caelo	Base cream, hydrophobic DAC	15221401	6
Caelo	Base cream, hydrophobic DAC	1522301	1
Caelo	Base cream, hydrophobic DAC	15240801	1
AHD/Pharmachem	Base cream, hydrophobic DAC	153625	4
Pharmachen GmbH/Gehe	Base cream, hydrophobic DAC	153625	2
Pharmachem	Base cream, hydrophobic DAC	153626	3
Pößneck/Gehe	Base cream, hydrophobic DAC	153626	1
AHD/Pharmachem	Base cream, hydrophobic DAC	153626	1
Pharmachem /Phönix	Base cream, hydrophobic DAC	153627	1
AHD/Pharmachem GmbH	Base cream, hydrophobic DAC	153627	1
Caelo	Base cream, hydrophobic DAC	153627	1
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	153627	1
Pharma/Phönix	Base cream, hydrophobic DAC	153627	1
Sanacorp	Base cream, hydrophobic DAC	15421201	3
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	15421201	1
Caelo	Base cream, hydrophobic DAC	15421303	1
Fagron	Base cream, hydrophobic DAC	15421401	1
Fagron	Base cream, hydrophobic DAC	15D03-T05-021168	1
Fagron	Base cream, hydrophobic DAC	15J13-T01-036887	1
Fagron	Base cream, hydrophobic DAC	15J13-T02-026891	1
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	1600825	1
pharmachem	Base cream, hydrophobic DAC	160130	1
Pharmachen GmbH/Gehe	Base cream, hydrophobic DAC	160130	2
Pharmachem GmbH & Co.KG	Base cream, hydrophobic DAC	160130	1
Pharma/Phönix	Base cream, hydrophobic DAC	160136	2
pharmachem	Base cream, hydrophobic DAC	160136	1
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	160136	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	160136	1
phönix	Base cream, hydrophobic DAC	160136	1
Bombastus	Base cream, hydrophobic DAC	160136	1
Pharmachem	Base cream, hydrophobic DAC	160137	1
Noweda/Pharmachem	Base cream, hydrophobic DAC	160137	1
Pharmachem	Base cream, hydrophobic DAC	160138	1
AHD / Pharmachem Pößneck	Base cream, hydrophobic DAC	160138	1
Caesar & Loretz/Noweda	Base cream, hydrophobic DAC	16074802	1
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	16074901	1
Caesar & Loretz/AHD	Base cream, hydrophobic DAC	16074901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz	Base cream, hydrophobic DAC	16075201	1
Pharmachem/ Gehe	Base cream, hydrophobic DAC	160814	1
AHD / Pharmachem Pößneck	Base cream, hydrophobic DAC	160814	1
Pharmachem/Phönix	Base cream, hydrophobic DAC	160814	1
Pharmachem GmbH & CoKG	Base cream, hydrophobic DAC	160814	1
Pharmachem	Base cream, hydrophobic DAC	160814	1
phönix	Base cream, hydrophobic DAC	160825	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	160825	2
Pharmachem GmbH Pößneck	Base cream, hydrophobic DAC	160825	1
Pharmachem Pößneck/Phönix	Base cream, hydrophobic DAC	160825	1
pharmachem	Base cream, hydrophobic DAC	160825	1
	Base cream, hydrophobic DAC	160825	1
Anzag /Pharmachem	Base cream, hydrophobic DAC	160825	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	160825	1
Pharmachem	Base cream, hydrophobic DAC	160826	3
Noweda	Base cream, hydrophobic DAC	16165403	1
Euro OTC	Base cream, hydrophobic DAC	16165501	1
caesar&Loretz GmbH,Phön. . .	Base cream, hydrophobic DAC	16165701	1
PHÖNIX	Base cream, hydrophobic DAC	16165701	1
Phoenix	Base cream, hydrophobic DAC	16174201	1
Pharmachen Pößneck	Base cream, hydrophobic DAC	162167	1
Pharmachem/Kehr	Base cream, hydrophobic DAC	162168	1
pharmachem	Base cream, hydrophobic DAC	162168	1
Pharmachem GmbH & Co. KG	Base cream, hydrophobic DAC	162168	2
Gehe-Pharmachem	Base cream, hydrophobic DAC	162168	1
Pharma/Phönix	Base cream, hydrophobic DAC	162168	1
Caelo	Base cream, hydrophobic DAC	162168	1
Pharmachen GmbH/Allianc. . .	Base cream, hydrophobic DAC	162169	1
Bombstus/San	Base cream, hydrophobic DAC	162168	1
Pharmachem/ Gehe	Base cream, hydrophobic DAC	162169	1
AHD / Pharmachem Pößneck	Base cream, hydrophobic DAC	162169	1
Pharmachem	Base cream, hydrophobic DAC	162170	4
Pößneck Pharma / Gehe	Base cream, hydrophobic DAC	162169	1
AHD/Pharmachem	Base cream, hydrophobic DAC	162170	1
pharmachem	Base cream, hydrophobic DAC	162170	1
Pharmachen	Base cream, hydrophobic DAC	162170	1
Sanacorp	Base cream, hydrophobic DAC	16293701	1
Sanacorp WE 25.02.17 -1. . .	Base cream, hydrophobic DAC	16293701	1
Caelo	Base cream, hydrophobic DAC	16293901	4
Noweda	Base cream, hydrophobic DAC	16293701	1
Phönix	Base cream, hydrophobic DAC	16293901	1
Pharmachem	Base cream, hydrophobic DAC	163237	2
Pharmachem/Phönix	Base cream, hydrophobic DAC	163238	1
Pharmachem	Base cream, hydrophobic DAC	163238	1
Caelo	Base cream, hydrophobic DAC	16330301	2
Sanacorp	Base cream, hydrophobic DAC	16330301	1
Caelo	Base cream, hydrophobic DAC	16330302	1
Fagron	Base cream, hydrophobic DAC	16B10-T01-030202	1
Fagron	Base cream, hydrophobic DAC	16E03T02032868	1
Fagron	Base cream, hydrophobic DAC	16E03-T02-032869	1
Fagron	Base cream, hydrophobic DAC	16E03T02-032868	1
Fagron	Base cream, hydrophobic DAC	16F27-T22-034684	3
Pößneck/Gehe	Base cream, hydrophobic DAC	170008	1
Pharmachem	Base cream, hydrophobic DAC	170008	1
Pharmachem/Phönix	Base cream, hydrophobic DAC	170008	1
Pharmachem / AHD	Base cream, hydrophobic DAC	170009	1
Caelo	Base cream, hydrophobic DAC	170009	1
Ebert/Pharmachem	Base cream, hydrophobic DAC	170008	1
Pharmachem GmbH & Co. KG	Base cream, hydrophobic DAC	170009	1
Allianc. Health. / Phar. . .	Base cream, hydrophobic DAC	170010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Pharmachem GmbH & Co. KG	Base cream, hydrophobic DAC	170010	1
Pharmachem	Base cream, hydrophobic DAC	170010	1
Pharmachem/ Phönix	Base cream, hydrophobic DAC	170010	2
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	170010	1
Alliance	Base cream, hydrophobic DAC	170010	2
Caelo	Base cream, hydrophobic DAC	17001093	1
Pharmachem GmbH	Base cream, hydrophobic DAC	17009	1
Caelo	Base cream, hydrophobic DAC		11
Caelo	Base cream, hydrophobic DAC	17073202	5
GEHE	Base cream, hydrophobic DAC	17073202	1
Noweda	Base cream, hydrophobic DAC	17073202	1
Caelo	Base cream, hydrophobic DAC	171117HB/17150401	1
Bombastus	Base cream, hydrophobic DAC	171126	2
Pharmachem	Base cream, hydrophobic DAC	171127	1
Bombastus	Base cream, hydrophobic DAC	171128	1
Pharmachem	Base cream, hydrophobic DAC	171128	2
Caesar & Loretz GmbH	Base cream, hydrophobic DAC	17112903	1
Caelo	Base cream, hydrophobic DAC	171219HB/17150401	1
Caelo	Base cream, hydrophobic DAC	17150401	3
Sanacorp W.E: 26.01.18 ...	Base cream, hydrophobic DAC	17150401	1
Gehe	Base cream, hydrophobic DAC	17150401	1
Caleo	Base cream, hydrophobic DAC	17150402	1
Caelo	Base cream, hydrophobic DAC	171815	1
Pharmachem GmbH/AHD	Base cream, hydrophobic DAC	171815	1
Caelo	Base cream, hydrophobic DAC	171816	1
Pharmachem	Base cream, hydrophobic DAC	171816	1
Pharmachem	Base cream, hydrophobic DAC	171817	1
Sanacorp/Pharmachem	Base cream, hydrophobic DAC	171817	1
Pharmachem/GEHE	Base cream, hydrophobic DAC	172194	4
Pharmachem/Noweda	Base cream, hydrophobic DAC	172194	1
Sanacorp/Pharmachem	Base cream, hydrophobic DAC	172195	2
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	172195	3
Pharmachem/Gehe	Base cream, hydrophobic DAC	172196	2
Pharmachen GmbH	Base cream, hydrophobic DAC	172195	1
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	172196	3
Pharmachem	Base cream, hydrophobic DAC	172196	1
Caelo	Base cream, hydrophobic DAC	17255301	1
Noweda	Base cream, hydrophobic DAC	17255301	2
Bombastus	Base cream, hydrophobic DAC	18002000	1
Caelo	Base cream, hydrophobic DAC	180308HB/17255301	1
Caelo	Base cream, hydrophobic DAC	18031505	1
GEHE	Base cream, hydrophobic DAC	18051101	1
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	180517	2
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	180755	1
Pharmachem/Sanacorp	Base cream, hydrophobic DAC	180757	7
Caelo	Base cream, hydrophobic DAC	18100801	1
phönix	Base cream, hydrophobic DAC	2/AKD12	1
Caelo	Base cream, hydrophobic DAC	20160222	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	231213	1
Caelo	Base cream, hydrophobic DAC	26011603	1
Bombastus	Base cream, hydrophobic DAC	284171	1
Bombastus	Base cream, hydrophobic DAC	285434	1
Bombastus	Base cream, hydrophobic DAC	285435	2
Bombastus	Base cream, hydrophobic DAC	285437	2
Bombastus	Base cream, hydrophobic DAC	288497	4
Bombastus	Base cream, hydrophobic DAC	289721	1
Caelo	Base cream, hydrophobic DAC	289720	1
Bombastus	Base cream, hydrophobic DAC	293162	3
Caelo	Base cream, hydrophobic DAC	293162	1
Bombastus	Base cream, hydrophobic DAC	293167	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Base cream, hydrophobic DAC	294685	6
Bombastus	Base cream, hydrophobic DAC	294686	3
Bombastus	Base cream, hydrophobic DAC	297046	5
Bombastus	Base cream, hydrophobic DAC	297044	1
Bombastus	Base cream, hydrophobic DAC	297047	6
Bombastus	Base cream, hydrophobic DAC	298351	1
Bombastus	Base cream, hydrophobic DAC	298354	2
Bombastus	Base cream, hydrophobic DAC	302133	1
Caelo	Base cream, hydrophobic DAC	302134	1
GEHE	Base cream, hydrophobic DAC	303373	1
Bombastus	Base cream, hydrophobic DAC	303374	1
Bombastus	Base cream, hydrophobic DAC	303377	2
Bombastus	Base cream, hydrophobic DAC	308509	1
GEHE	Base cream, hydrophobic DAC	310824	1
Bombastus	Base cream, hydrophobic DAC	310823	1
Bombastus	Base cream, hydrophobic DAC	310826	1
Caelo	Base cream, hydrophobic DAC	5013q-03497	1
Bombastus	Base cream, hydrophobic DAC	60080401	2
Caelo	Base cream, hydrophobic DAC	5421201	1
Sanacorp	Base cream, hydrophobic DAC	4/AKD3	1
Caelo	Base cream, hydrophobic DAC	46074802	1
Caelo	Base cream, hydrophobic DAC	9AKD15	1
Bombastus	Base cream, hydrophobic DAC		1
Bombastus	Base cream, hydrophobic DAC	269276	1
Bombastus	Base cream, hydrophobic DAC	271091	1
Bombastus	Base cream, hydrophobic DAC	271093	1
Ahd/PharmachemPößneck	Base cream, hydrophobic DAC	21013	1
Bombastus	Base cream, hydrophobic DAC	271095	1
Pharmachem/30	Base cream, hydrophobic DAC	21013	1
Pharmachen	Base cream, hydrophobic DAC	21013	1
Fagron	Base cream, hydrophobic DAC	513A-03497	1
Caelo	Base cream, hydrophobic DAC	2111201	1
Phoenix/Pharmachem	Base cream, hydrophobic DAC	110413	1
PharmaChem	Base cream, hydrophobic DAC	110413	1
pharmachem	Base cream, hydrophobic DAC	110413	1
Caelo	Base cream, hydrophobic DAC	12064047	1
Pharmachem/Gehe	Base cream, hydrophobic DAC	120712	2
PharmaChem / Noweda	Base cream, hydrophobic DAC	120712	1
Erika apotheke	Base cream, hydrophobic DAC	12151601	1
Noweda/ Ichtyol	Base cream, hydrophobic DAC	12151601	2
Sanacorp	Base cream, hydrophobic DAC	12151602	1
Caelo	Base cream, hydrophobic DAC	1215601	1
Caelo	Base cream, hydrophobic DAC	12217603	6
Noweda	Base cream, hydrophobic DAC	12229302	4
Bombastus	Base cream, hydrophobic DAC	12229302	1
Caelo	Base cream, hydrophobic DAC	12229304	2
Sanacorp	Base cream, hydrophobic DAC	12229304	1
Sanacorp	Base cream, hydrophobic DAC	12308102	2
AHD, 02.04.2013	Base cream, hydrophobic DAC	12308102	1
Caelo	Base cream, hydrophobic DAC	12308103	2
Caelo	Base cream, hydrophobic DAC	123082011	1
Phönix	Base cream, hydrophobic DAC	1230822	1
Fagron	Base cream, hydrophobic DAC	12D10-T02	1
Fagron	Base cream, hydrophobic DAC	12G026-0001	1
Kehr	Base cream, hydrophobic DAC	12Go9-T03	1
Fagron	Base cream, hydrophobic DAC	12G09-T03	2
Kehr	Base cream, hydrophobic DAC	12Gog-To3	1
Caesar & Loretz GmbH, 4. .	Base cream, hydrophobic DAC	13019303	1
Taoasis	Base cream, hydrophobic DAC	130719-3198	1
Caelo	Base cream, hydrophobic DAC	13085504	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Base cream, hydrophobic DAC	13085701	2
Gehe	Base cream, hydrophobic DAC	13085701	1
Caelo	Base cream, hydrophobic DAC	13085702	1
Caelo	Base cream, hydrophobic DAC	13086701	1
Caelo	Base cream, hydrophobic DAC	1309302	1
Pharmachem	Base cream, hydrophobic DAC	131211	1
Caelo	Base cream, hydrophobic DAC	13127502	2
Phoenix	Base cream, hydrophobic DAC	13175902	1
cealo	Base cream, hydrophobic DAC	13274902	1
Caelo	Base cream, hydrophobic DAC	13275302	7
Caesar&loretz,Phoenix	Base cream, hydrophobic DAC	13386401	1
Fagron	Base cream, hydrophobic DAC	13B04-T17-001040	2
Fagron	Base cream, hydrophobic DAC	13E13T01003619	1
Fagron	Base cream, hydrophobic DAC	13E13-T01-003620	3
Fagron	Base cream, hydrophobic DAC	13H13T06	1
Caelo	Base cream, hydrophobic DAC	1400601	1
pharmachem	Base cream, hydrophobic DAC	140886	1
Caesar&loretz,Phoenix	Base cream, hydrophobic DAC	15001002	1
Pharmachem	Base cream, hydrophobic DAC	141489	1
Caelo	Base cream, hydrophobic DAC	21122013B	1
Pharmachem.Pöfneck/Phoe...	Base cream, hydrophobic DAC	230713	1
Pharmachem / Anzag	Base cream, hydrophobic DAC	230713	2
Pharmachem	Base cream, hydrophobic DAC	230713	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	230713	1
PharmaChin/ Anzag	Base cream, hydrophobic DAC	230713	1
Pharmachen	Base cream, hydrophobic DAC	231213	1
pharmachem	Base cream, hydrophobic DAC	231213	1
Pharmachem/Noweda	Base cream, hydrophobic DAC	231213	1
Pharmachem / Anzag	Base cream, hydrophobic DAC	240113	1
Pharmachem GmbH	Base cream, hydrophobic DAC	261112	1
pharmachen	Base cream, hydrophobic DAC	240113	1
Pharmachen	Base cream, hydrophobic DAC	240113	1
Hersteller: Pharmachem/...	Base cream, hydrophobic DAC	261112	3
Pharmachem / Anzag	Base cream, hydrophobic DAC	261112	1
Pharmachen Pöfneck	Base cream, hydrophobic DAC	261112	1
Bombastus	Base cream, hydrophobic DAC	276087	1
Pharmachem Pöfneck 004	Base cream, hydrophobic DAC	261112	1
Gehe	Base cream, hydrophobic DAC	261112	1
Phönix Gotha	Base cream, hydrophobic DAC	277696	1
Bombastus	Base cream, hydrophobic DAC	277696	1
Bombastus	Base cream, hydrophobic DAC	279145	2
Bombastus	Base cream, hydrophobic DAC	279148	3
Bombastus	Base cream, hydrophobic DAC	282161	1
Bombastus	Base cream, hydrophobic DAC	280809	1
Bombastus	Base cream, hydrophobic DAC	284170	1
Caelo	Base cream, hydrophobic DAC	410211-1-P	1
Noweda	Base cream, hydrophobic DAC	11289203	1
Pöfneck Pharmachem	Base cream, hydrophobic DAC	131011	1
Bombastus	Hydrophilic cream, non-ionic...	2293552	1
GEHE/Apomix/25.03.2014	Hydrophilic cream, non-ionic...	1A140114	1
apomix/Anzag	Hydrophilic cream, non-ionic...	1a200116	1
apomix/Noweda	Hydrophilic cream, non-ionic...	1A210115	1
Apomix	Hydrophilic cream, non-ionic...	1A210115	1
gehe	Hydrophilic cream, non-ionic...	1B070114	1
Noweda	Hydrophilic cream, non-ionic...	1B070116	1
Apomix	Hydrophilic cream, non-ionic...	1B120115	1
Apomix/Noweda	Hydrophilic cream, non-ionic...	2A090215	1
apomix	Hydrophilic cream, non-ionic...	2A090215	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	2A090215	1
Gehe/Apomix	Hydrophilic cream, non-ionic...	2A140114	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hydrophilic cream, non-ionic...	2A140114	2
apomix/Gehe	Hydrophilic cream, non-ionic...	2B300117	1
Apomix/Anzag	Hydrophilic cream, non-ionic...	3A090215	1
AMH Niemann GmbH&Co KG	Hydrophilic cream, non-ionic...	2B130114	1
apomix/Phönix	Hydrophilic cream, non-ionic...	3A090215	1
Apomix/Gehe	Hydrophilic cream, non-ionic...	3B020215	1
Apomix/ Noweda	Hydrophilic cream, non-ionic...	3B020215	1
Apomix / Noweda	Hydrophilic cream, non-ionic...	3B020215	1
Apomix / Phönix	Hydrophilic cream, non-ionic...	4A250214	2
Apomix/Noweda	Hydrophilic cream, non-ionic...	4A240215	1
APOMIX/Noweda	Hydrophilic cream, non-ionic...	4a230216	1
apomix/Anzag	Hydrophilic cream, non-ionic...	4A240215	1
Gehe/Apomix/11.11.14	Hydrophilic cream, non-ionic...	4A250214	1
Raben-Apotheke Mildenau	Hydrophilic cream, non-ionic...	4B100214	1
Apomix	Hydrophilic cream, non-ionic...	4B170216	1
Apomix	Hydrophilic cream, non-ionic...	4B210217	1
Apomix	Hydrophilic cream, non-ionic...	5B01816	1
Apomix / Noweda	Hydrophilic cream, non-ionic...	5B150317	1
apomix/Gehe	Hydrophilic cream, non-ionic...	5B150317	1
Apomix/ Gehe	Hydrophilic cream, non-ionic...	5B290216	1
Apomix	Hydrophilic cream, non-ionic...	6B020315	1
Apomix / Sanacorp	Hydrophilic cream, non-ionic...	6B020315	1
Apomix	Hydrophilic cream, non-ionic...	6B100314	1
Pharmachem	Hydrophilic cream, non-ionic...	6B150316	1
Apomix/Noweda	Hydrophilic cream, non-ionic...	7A020414	1
Apomix	Hydrophilic cream, non-ionic...	7A020414	1
Apomix / Phönix	Hydrophilic cream, non-ionic...	7A020414	1
Gehe/apomix/09.10.14	Hydrophilic cream, non-ionic...	7A020414	1
Apomix	Hydrophilic cream, non-ionic...	7A140217	1
Apomix/Kehr	Hydrophilic cream, non-ionic...	7A140217	1
Allianc. Health. /Apomix	Hydrophilic cream, non-ionic...	7A310316	1
PKH GmbH/Noweda	Hydrophilic cream, non-ionic...	7B030417	1
PKH GmbH	Hydrophilic cream, non-ionic...	7B030417	1
apomix/Noweda	Hydrophilic cream, non-ionic...	7A310316	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	7a310316	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	8A120416	1
Apomix	Hydrophilic cream, non-ionic...	8A120416	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	8B060416	1
Apomix/ Sanacorp	Hydrophilic cream, non-ionic...	8B030417	1
Apomix	Hydrophilic cream, non-ionic...	8B070414	1
apomix/Phönix	Hydrophilic cream, non-ionic...	8B070414	1
Pharmachem/Phoenix	Hydrophilic cream, non-ionic...	90714	1
Alliance	Hydrophilic cream, non-ionic...	90714	1
Noweda	Hydrophilic cream, non-ionic...	90714	1
Apomix/ Sanacorp	Hydrophilic cream, non-ionic...	9A210317	1
apomix/Anzag	Hydrophilic cream, non-ionic...	9a210317	1
Apomix	Hydrophilic cream, non-ionic...	9B020516	1
Apomix / Noweda	Hydrophilic cream, non-ionic...	9B200415	1
apomix	Hydrophilic cream, non-ionic...	10B190516	1
Apomix	Hydrophilic cream, non-ionic...	10B190516	1
PKH GmbH/ NOWEDA	Hydrophilic cream, non-ionic...	10B200415	1
Phönix Gotha	Hydrophilic cream, non-ionic...	10B200415	1
Apomix	Hydrophilic cream, non-ionic...	11B110515	1
Apomix/Noweda	Hydrophilic cream, non-ionic...	10B200415	1
PKH Halle	Hydrophilic cream, non-ionic...	11B310517	1
Apomix / Noweda	Hydrophilic cream, non-ionic...	11B310517	1
Fagron	Hydrophilic cream, non-ionic...	12A020315	1
Fagron	Hydrophilic cream, non-ionic...	12A180416	2
PKH GmbH	Hydrophilic cream, non-ionic...	12B1306161	1
Bombastus	Hydrophilic cream, non-ionic...	12B130626	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix / Sanacorp	Hydrophilic cream, non-ionic...	12B210514	1
Raben-Apotheke Mildenau	Hydrophilic cream, non-ionic...	12B210514	1
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	133911	2
Leyh Pharma	Hydrophilic cream, non-ionic...	134841	1
Fagron	Hydrophilic cream, non-ionic...	13A180416	1
Apomix	Hydrophilic cream, non-ionic...	13B080615	2
apomix	Hydrophilic cream, non-ionic...	13B080615	1
apomix	Hydrophilic cream, non-ionic...	13B280616	2
Apomix / Noweda	Hydrophilic cream, non-ionic...	13B280617	2
Fagron	Hydrophilic cream, non-ionic...	13K15-M09	1
Fagron	Hydrophilic cream, non-ionic...	13K20-M09	1
Fagron	Hydrophilic cream, non-ionic...	13K20-M10	1
Fagron	Hydrophilic cream, non-ionic...	13K20-M11	1
Leyh Pharma/Gehe	Hydrophilic cream, non-ionic...	140551	1
Caelo	Hydrophilic cream, non-ionic...	14063003	1
Caelo	Hydrophilic cream, non-ionic...	14137002	5
Caelo	Hydrophilic cream, non-ionic...	14221101	2
Phönix/Leyh Pharma	Hydrophilic cream, non-ionic...	142311	1
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	142311	1
Leyh Pharma / Phönix	Hydrophilic cream, non-ionic...	142311	1
Leyh Pharma/Phönix	Hydrophilic cream, non-ionic...	143011	1
Leyh/Gehe	Hydrophilic cream, non-ionic...	143011	2
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	143011	1
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	144021	1
Fagron	Hydrophilic cream, non-ionic...	14A190416	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	14A220517	1
Apomix/Kehr	Hydrophilic cream, non-ionic...	14B120716	1
Fagron	Hydrophilic cream, non-ionic...	14d28-m01	1
Fagron	Hydrophilic cream, non-ionic...	14D28-M02	1
Fagron	Hydrophilic cream, non-ionic...	14D28-M03	1
Fagron	Hydrophilic cream, non-ionic...	14D28-M05	1
Fagron	Hydrophilic cream, non-ionic...	14D28-M05	2
Fagron	Hydrophilic cream, non-ionic...	14D28-Mo9	1
Fagron	Hydrophilic cream, non-ionic...	14E07-M10	1
Fagron	Hydrophilic cream, non-ionic...	14E14-M06	1
Fagron	Hydrophilic cream, non-ionic...	14E07	1
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	150341	1
Phönix Gotha	Hydrophilic cream, non-ionic...	150341	1
Gehe/Leyh	Hydrophilic cream, non-ionic...	150341	1
Leyh Pharma/ Phönix	Hydrophilic cream, non-ionic...	150751	1
Caelo	Hydrophilic cream, non-ionic...	15325802	1
Caelo	Hydrophilic cream, non-ionic...	15365601	2
Fagron	Hydrophilic cream, non-ionic...	15A-200416	1
Apomix/ Phönix	Hydrophilic cream, non-ionic...	15B070714	1
Fagron	Hydrophilic cream, non-ionic...	15A200416	1
AMH-NiemannGmbH&CO,15.0. ...	Hydrophilic cream, non-ionic...	15B010817	1
Fagron	Hydrophilic cream, non-ionic...	15C04-M06	1
PKH Halle	Hydrophilic cream, non-ionic...	15B080816	1
Apomix/ Phönix	Hydrophilic cream, non-ionic...	15B090713	1
Bombastus	Hydrophilic cream, non-ionic...	160824UEASR	1
Caelo	Hydrophilic cream, non-ionic...	16169501	3
SanacorpWE:01.12.2016 ...	Hydrophilic cream, non-ionic...	16222001	1
Caelo	Hydrophilic cream, non-ionic...	16222102	2
Leyh Pharma/Noweda	Hydrophilic cream, non-ionic...	162241	1
Noweda	Hydrophilic cream, non-ionic...	16321701	1
Leyh Pharma/Noweda	Hydrophilic cream, non-ionic...	163611	1
Gehe	Hydrophilic cream, non-ionic...	17000737	1
GEHE	Hydrophilic cream, non-ionic...	17000738	1
Bombastus	Hydrophilic cream, non-ionic...	17000741	1
Bombastus	Hydrophilic cream, non-ionic...	17000745	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrophilic cream, non-ionic...	17000746	1
Bombastus	Hydrophilic cream, non-ionic...	17000749	1
Bombastus	Hydrophilic cream, non-ionic...	17000750	2
Bombastus	Hydrophilic cream, non-ionic...	17002159	1
Bombastus	Hydrophilic cream, non-ionic...	1700750	1
Caelo	Hydrophilic cream, non-ionic...	170427NHCSR	1
Leyh Pharma/Noweda	Hydrophilic cream, non-ionic...	170951	1
Caelo	Hydrophilic cream, non-ionic...	17144001	2
Caelo	Hydrophilic cream, non-ionic...	17212601	3
Caelo	Hydrophilic cream, non-ionic...	17294002	2
Caelo	Hydrophilic cream, non-ionic...	1744001	1
Apomix	Hydrophilic cream, non-ionic...	17B120814	2
Caelo	Hydrophilic cream, non-ionic...	18001330	1
Bombastus	Hydrophilic cream, non-ionic...	18001332	1
Bombastus	Hydrophilic cream, non-ionic...	18001334	1
Bombastus	Hydrophilic cream, non-ionic...	18002063	1
Caelo	Hydrophilic cream, non-ionic...	180306UEN/17313902	1
Gehe/Apomix	Hydrophilic cream, non-ionic...	18A190315	1
Raben-Apotheke Mildenau	Hydrophilic cream, non-ionic...	18B020914	2
Apomix / Sanacorp	Hydrophilic cream, non-ionic...	18B240815	1
apomix noweder 23.10.2015	Hydrophilic cream, non-ionic...	18B240815	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	19A300315	2
apomix	Hydrophilic cream, non-ionic...	19A300315	1
apomix/Phönix	Hydrophilic cream, non-ionic...	19A300315	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	19A300315	1
Apomix/NOWEDA	Hydrophilic cream, non-ionic...	19A300315	1
Sanacorp/ apomix	Hydrophilic cream, non-ionic...	19B021017	1
Caelo	Hydrophilic cream, non-ionic...	201015	2
apomix	Hydrophilic cream, non-ionic...	20A280415	1
PKH Halle/Noweda	Hydrophilic cream, non-ionic...	20A280415	1
Apomix/Gehe	Hydrophilic cream, non-ionic...	20B210915	1
Apomix/ Phönix	Hydrophilic cream, non-ionic...	20B230913	1
Rabenapotheker	Hydrophilic cream, non-ionic...	20B290914	2
Phönix Gotha	Hydrophilic cream, non-ionic...	21B141014	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	21b141116	1
PKH GmbH	Hydrophilic cream, non-ionic...	21B141116	2
Caelo	Hydrophilic cream, non-ionic...	21B141116	1
apomix	Hydrophilic cream, non-ionic...	21B14116	1
Pharmachem/Phoenix	Hydrophilic cream, non-ionic...	220713	12
Pharmachem / Phönix	Hydrophilic cream, non-ionic...	220713	1
apomix/Phönix	Hydrophilic cream, non-ionic...	22B061114	1
apomix/Phönix	Hydrophilic cream, non-ionic...	23A100615	1
Apomix /Alliance	Hydrophilic cream, non-ionic...	22B281013	1
Apomix	Hydrophilic cream, non-ionic...	23A280416	2
apomix	Hydrophilic cream, non-ionic...	24B241115	1
Pharmachem/Noweda	Hydrophilic cream, non-ionic...	240217	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	26A080715	1
Apomix/Anzag	Hydrophilic cream, non-ionic...	26A080715	1
Apomix/Sana	Hydrophilic cream, non-ionic...	25B081215	1
apomix / Anzag	Hydrophilic cream, non-ionic...	26A210616	1
Bombastus	Hydrophilic cream, non-ionic...	276593	1
PKH/Noweda	Hydrophilic cream, non-ionic...	27A220715	1
gehe	Hydrophilic cream, non-ionic...	280990	1
Bombastus	Hydrophilic cream, non-ionic...	282212	3
Bombastus	Hydrophilic cream, non-ionic...	282215	2
Caesar & Loretz GmbH Hi...	Hydrophilic cream, non-ionic...	282216	1
Caelo	Hydrophilic cream, non-ionic...	283634	1
Bombastus	Hydrophilic cream, non-ionic...	283634	1
Bombastus	Hydrophilic cream, non-ionic...	283635	2
Gehe	Hydrophilic cream, non-ionic...	283637	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
gehe	Hydrophilic cream, non-ionic...	283693	1
Bombastus	Hydrophilic cream, non-ionic...	283693	1
Bombastus	Hydrophilic cream, non-ionic...	283694	4
Bombastus	Hydrophilic cream, non-ionic...	284891	3
Bombastus	Hydrophilic cream, non-ionic...	284893	2
Bombastus	Hydrophilic cream, non-ionic...	284895	5
Bombastus	Hydrophilic cream, non-ionic...	284896	2
Bombastus	Hydrophilic cream, non-ionic...	284897	1
Bombastus	Hydrophilic cream, non-ionic...	284898	2
Herbasinica	Hydrophilic cream, non-ionic...	284899	1
Nachprüfung	Hydrophilic cream, non-ionic...	284899	2
Bombastus	Hydrophilic cream, non-ionic...	286323	1
Hedinger	Hydrophilic cream, non-ionic...	286326	1
Caelo	Hydrophilic cream, non-ionic...	286327	1
Bombastus	Hydrophilic cream, non-ionic...	286328	3
Bombastus	Hydrophilic cream, non-ionic...	286329	2
Bombastus	Hydrophilic cream, non-ionic...	286838	2
Bombastus	Hydrophilic cream, non-ionic...	287710	2
Gehe 02.09.14	Hydrophilic cream, non-ionic...	287709	1
Bombastus	Hydrophilic cream, non-ionic...	287711	9
Bombastus	Hydrophilic cream, non-ionic...	287712	2
AHD 08.09.14 EK 5,31	Hydrophilic cream, non-ionic...	287711	1
Holdermann	Hydrophilic cream, non-ionic...	287712	1
Bombastus	Hydrophilic cream, non-ionic...	287713	6
Bombastus	Hydrophilic cream, non-ionic...	287717	2
Bombastus	Hydrophilic cream, non-ionic...	287762	2
Bombastus	Hydrophilic cream, non-ionic...	288793	1
Bombastus	Hydrophilic cream, non-ionic...	288796	1
Bombastus	Hydrophilic cream, non-ionic...	288797	1
Bombastus	Hydrophilic cream, non-ionic...	288882	3
Bombastus	Hydrophilic cream, non-ionic...	288798	1
Bombastus	Hydrophilic cream, non-ionic...	289326	1
Bombastus	Hydrophilic cream, non-ionic...	289937	1
Gehe	Hydrophilic cream, non-ionic...	289941	1
Bombastus	Hydrophilic cream, non-ionic...	289940	2
Gehe	Hydrophilic cream, non-ionic...	289937	1
Bombastus	Hydrophilic cream, non-ionic...	289941	1
Bombastus	Hydrophilic cream, non-ionic...	289942	5
Bombastus	Hydrophilic cream, non-ionic...	289943	4
Caelo	Hydrophilic cream, non-ionic...	289944	1
apomix	Hydrophilic cream, non-ionic...	28A060815	1
Bombastus	Hydrophilic cream, non-ionic...	291382	1
Noweda	Hydrophilic cream, non-ionic...	291385	1
Euro OTC	Hydrophilic cream, non-ionic...	291388	1
Bombastus	Hydrophilic cream, non-ionic...	291388	9
Caelo	Hydrophilic cream, non-ionic...	291388	1
Bombastus	Hydrophilic cream, non-ionic...	292737	1
Bombastus	Hydrophilic cream, non-ionic...	292738	1
Euro OTC	Hydrophilic cream, non-ionic...	292737	1
Bombastus	Hydrophilic cream, non-ionic...	292740	1
Noweda	Hydrophilic cream, non-ionic...	292805	1
Bombastus	Hydrophilic cream, non-ionic...	292806	2
Bombastus	Hydrophilic cream, non-ionic...	292808	2
Bombastus	Hydrophilic cream, non-ionic...	292809	2
Bombastus	Hydrophilic cream, non-ionic...	292810	1
Anzag	Hydrophilic cream, non-ionic...	294311	1
Bombastus	Hydrophilic cream, non-ionic...	294305	1
Bombastus	Hydrophilic cream, non-ionic...	294311	3
Caelo	Hydrophilic cream, non-ionic...	294310	1
Noweda	Hydrophilic cream, non-ionic...	294311	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrophilic cream, non-ionic...	294340	1
Bombastus	Hydrophilic cream, non-ionic...	294313	1
Bombastus	Hydrophilic cream, non-ionic...	294390	1
Noweda	Hydrophilic cream, non-ionic...	294390	2
Bombastus	Hydrophilic cream, non-ionic...	294408	3
Beiersdorf/Kehr	Hydrophilic cream, non-ionic...	294407	1
Noweda	Hydrophilic cream, non-ionic...	294408	1
Bombastus	Hydrophilic cream, non-ionic...	295644	1
Bombastus / Noweda	Hydrophilic cream, non-ionic...	295647	1
Bombastus	Hydrophilic cream, non-ionic...	295591	1
Bombastus	Hydrophilic cream, non-ionic...	295648	3
Noweda	Hydrophilic cream, non-ionic...	295648	1
Bombastus	Hydrophilic cream, non-ionic...	295649	3
Caelo	Hydrophilic cream, non-ionic...	295649	1
Bombastus	Hydrophilic cream, non-ionic...	295651	1
Bombastus	Hydrophilic cream, non-ionic...	295650	1
Gehe	Hydrophilic cream, non-ionic...	295653	1
Bombastus	Hydrophilic cream, non-ionic...	295653	1
Bombastus	Hydrophilic cream, non-ionic...	295654	1
Bombastus	Hydrophilic cream, non-ionic...	295655	8
Bombastus	Hydrophilic cream, non-ionic...	297299	1
Bombastus	Hydrophilic cream, non-ionic...	297300	1
Bombastus	Hydrophilic cream, non-ionic...	297303	2
Bombastus	Hydrophilic cream, non-ionic...	297308	10
Noweda	Hydrophilic cream, non-ionic...	297309	1
Bombastus	Hydrophilic cream, non-ionic...	297371	3
Bombastus	Hydrophilic cream, non-ionic...	29739	1
Caelo	Hydrophilic cream, non-ionic...	297374	2
Bombastus	Hydrophilic cream, non-ionic...	298979	1
Bombastus	Hydrophilic cream, non-ionic...	298982	2
Bombastus	Hydrophilic cream, non-ionic...	298983	4
Bombastus	Hydrophilic cream, non-ionic...	298984	2
apomix / Anzag	Hydrophilic cream, non-ionic...	29A020816	2
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	29A110815	1
Bombastus	Hydrophilic cream, non-ionic...	2x295649	1
Bombastus	Hydrophilic cream, non-ionic...	300035	4
Phoenix	Hydrophilic cream, non-ionic...	300035	1
Bombastus	Hydrophilic cream, non-ionic...	300759	3
Bombastus	Hydrophilic cream, non-ionic...	300762	1
Noweda	Hydrophilic cream, non-ionic...	300760	1
Bombastus	Hydrophilic cream, non-ionic...	300760	1
Bombastus	Hydrophilic cream, non-ionic...	300765	1
Bombastus	Hydrophilic cream, non-ionic...	30162	1
Bombastus	Hydrophilic cream, non-ionic...	301621	2
Bombastus	Hydrophilic cream, non-ionic...	301626	3
Caelo	Hydrophilic cream, non-ionic...	301627	1
Bombastus	Hydrophilic cream, non-ionic...	302683	1
Bombastus	Hydrophilic cream, non-ionic...	302685	2
GEHE	Hydrophilic cream, non-ionic...	302685	1
Bombastus	Hydrophilic cream, non-ionic...	302687	1
Bombastus	Hydrophilic cream, non-ionic...	302689	1
Bombastus	Hydrophilic cream, non-ionic...	302691	3
Bombastus	Hydrophilic cream, non-ionic...	303764	3
Bombastus	Hydrophilic cream, non-ionic...	304778	3
Bombastus	Hydrophilic cream, non-ionic...	304779	1
Bombastus	Hydrophilic cream, non-ionic...	305143	3
Bombastus	Hydrophilic cream, non-ionic...	306124	5
GEHE	Hydrophilic cream, non-ionic...	306124	1
Caelo	Hydrophilic cream, non-ionic...	306128	1
Bombastus	Hydrophilic cream, non-ionic...	306128	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrophilic cream, non-ionic...	306129	1
Bombastus	Hydrophilic cream, non-ionic...	30683	1
GEHE	Hydrophilic cream, non-ionic...	307521	1
Bombastus	Hydrophilic cream, non-ionic...	307523	1
Bombastus	Hydrophilic cream, non-ionic...	307525	2
Bombastus	Hydrophilic cream, non-ionic...	307524	1
Bombastus	Hydrophilic cream, non-ionic...	307528	1
Bombastus	Hydrophilic cream, non-ionic...	307529	1
Caelo	Hydrophilic cream, non-ionic...	309903	1
Ebert-Jacobi	Hydrophilic cream, non-ionic...	309904	1
Bombastus	Hydrophilic cream, non-ionic...	309905	1
Bombastus	Hydrophilic cream, non-ionic...	309907	2
GEHE/Apomix/04.02.14	Hydrophilic cream, non-ionic...	30A030613	2
Bombastus	Hydrophilic cream, non-ionic...	310514	1
Bombastus	Hydrophilic cream, non-ionic...	310516	1
Bombastus	Hydrophilic cream, non-ionic...	310517	2
Bombastus	Hydrophilic cream, non-ionic...	311132	1
Bombastus	Hydrophilic cream, non-ionic...	311133	1
Apomix/Kehr	Hydrophilic cream, non-ionic...	32A210916	1
Apomix/AHD/24.06.14	Hydrophilic cream, non-ionic...	32A250414	1
apomix	Hydrophilic cream, non-ionic...	32A220915	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	32A220915	1
apomix/Holdermann	Hydrophilic cream, non-ionic...	32A250414	2
Fagron	Hydrophilic cream, non-ionic...	33A180912	2
AMH Niemann GmbH & Co. KG	Hydrophilic cream, non-ionic...	34A210916	1
apomix / Anzag	Hydrophilic cream, non-ionic...	35A251016	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	36A140514	1
AHD/apomix/23.07.14	Hydrophilic cream, non-ionic...	36A140514	1
Allianc. Health. / Apomix	Hydrophilic cream, non-ionic...	36A261016	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	37A071216	1
Apomix	Hydrophilic cream, non-ionic...	36A140514a	1
Apomix /Noweda	Hydrophilic cream, non-ionic...	36A261016	1
Raben-Apotheke Mildenau	Hydrophilic cream, non-ionic...	37A270514	1
PKH GmbH	Hydrophilic cream, non-ionic...	36A261016	1
Apomix	Hydrophilic cream, non-ionic...	38A050614	1
Alliance Healthcare/Apo...	Hydrophilic cream, non-ionic...	38A050614	1
Fagron	Hydrophilic cream, non-ionic...	38A121216	1
Apomix/Noweda	Hydrophilic cream, non-ionic...	39A160614	1
Apomix	Hydrophilic cream, non-ionic...	40A160614	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	40A160614	1
Fagron	Hydrophilic cream, non-ionic...	41A141216	1
Anzag/apomix/26.03.14	Hydrophilic cream, non-ionic...	40A281013	1
Alliance Healthcare/Apo...	Hydrophilic cream, non-ionic...	42A070714	1
Fagron	Hydrophilic cream, non-ionic...	4315A-03195	1
apomix/Holdermann	Hydrophilic cream, non-ionic...	43A140714	1
Fagron	Hydrophilic cream, non-ionic...	43A191216	1
Fagron	Hydrophilic cream, non-ionic...	46A061113	1
Frago/Noweda	Hydrophilic cream, non-ionic...	46A201216	1
Apomix/Anzag	Hydrophilic cream, non-ionic...	47A281015	1
apomix/Phönix	Hydrophilic cream, non-ionic...	47A281015	1
Fagron	Hydrophilic cream, non-ionic...	48A071113	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	47A281015	1
Apomix	Hydrophilic cream, non-ionic...	48A081014	1
Gehe/apomix/11.12.14	Hydrophilic cream, non-ionic...	48A081014	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	48A081014	1
Apomix/Kehr	Hydrophilic cream, non-ionic...	48A181115	1
Apomix	Hydrophilic cream, non-ionic...	49A241115	1
Alliance Healthcare/Apo...	Hydrophilic cream, non-ionic...	50A101114	1
Apomix	Hydrophilic cream, non-ionic...	50A101114	1
Caelo	Hydrophilic cream, non-ionic...	50A111113	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix/Noweda	Hydrophilic cream, non-ionic...	50A161215	1
Fagron	Hydrophilic cream, non-ionic...	51A121113	2
Apomix	Hydrophilic cream, non-ionic...	51A261114	1
apomix/Gehe	Hydrophilic cream, non-ionic...	51A261114	1
Fagron	Hydrophilic cream, non-ionic...	5214I-03195	1
apomix	Hydrophilic cream, non-ionic...	52a081214	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	52A081214	1
Apomix / AHD	Hydrophilic cream, non-ionic...	53A021213	1
Apomix/Sanacorp	Hydrophilic cream, non-ionic...	53A191113	1
Apomix/Noweda	Hydrophilic cream, non-ionic...	54A021213	1
Phönix	Hydrophilic cream, non-ionic...		1
Bombastus	Hydrophilic cream, non-ionic...	253825	1
Bombastus	Hydrophilic cream, non-ionic...	264427	1
Bombastus	Hydrophilic cream, non-ionic...	267068	1
Bombastus	Hydrophilic cream, non-ionic...	267185	2
Caelo	Hydrophilic cream, non-ionic...	267186	1
Bombastus	Hydrophilic cream, non-ionic...	268349	2
Bombastus	Hydrophilic cream, non-ionic...	268353	1
Bombastus	Hydrophilic cream, non-ionic...	269712	1
Bombastus	Hydrophilic cream, non-ionic...	269711	1
Bombastus	Hydrophilic cream, non-ionic...	269716	1
Bombastus	Hydrophilic cream, non-ionic...	269715	1
Bombastus	Hydrophilic cream, non-ionic...	271473	1
Bombastus	Hydrophilic cream, non-ionic...	271518	1
Gehe	Hydrophilic cream, non-ionic...	271520	2
Bombastus	Hydrophilic cream, non-ionic...	271520	1
Bombastus	Hydrophilic cream, non-ionic...	271521	2
Gehe	Hydrophilic cream, non-ionic...	271521	1
Bombastus	Hydrophilic cream, non-ionic...	272472	1
Bombastus	Hydrophilic cream, non-ionic...	273132	1
Bombastus	Hydrophilic cream, non-ionic...	273137	1
Bombastus	Hydrophilic cream, non-ionic...	274212	2
Bombastus	Hydrophilic cream, non-ionic...	274213	1
Bombastus	Hydrophilic cream, non-ionic...	274214	1
Gehe	Hydrophilic cream, non-ionic...	274215	1
Bombastus	Hydrophilic cream, non-ionic...	274217	2
Bombastus	Hydrophilic cream, non-ionic...	274222	1
Apomix/Alliance Healthc...	Hydrophilic cream, non-ionic...	8/14	1
Apomix	Hydrophilic cream, non-ionic...	1A030112	1
Apomix	Hydrophilic cream, non-ionic...	1A030113	1
Alliance Healthcare	Hydrophilic cream, non-ionic...	1A030113	2
apomix	Hydrophilic cream, non-ionic...	2A140113	1
APOMIX /11	Hydrophilic cream, non-ionic...	1B070114	1
apomix/GEHE/22.05.2013	Hydrophilic cream, non-ionic...	2A140113	1
Apomix	Hydrophilic cream, non-ionic...	2A140113	1
apomix	Hydrophilic cream, non-ionic...	2B140113	1
Apomix	Hydrophilic cream, non-ionic...	2B140113	2
Phönix/ Apomix	Hydrophilic cream, non-ionic...	2B140113	1
Phönix/ Apomix	Hydrophilic cream, non-ionic...	2B160412	1
apomix	Hydrophilic cream, non-ionic...	3A240413	1
Phönix/ Apomix	Hydrophilic cream, non-ionic...	3B230112	1
Pharmachem/Phönix	Hydrophilic cream, non-ionic...	50912	4
apo	Hydrophilic cream, non-ionic...	4807013	1
Apomix	Hydrophilic cream, non-ionic...	4B070213	1
Fagron	Hydrophilic cream, non-ionic...	5A230113	1
Gehe	Hydrophilic cream, non-ionic...	5B250213	1
Apomix	Hydrophilic cream, non-ionic...	5B250213	1
Apomix	Hydrophilic cream, non-ionic...	6B110313	2
AMH Niemann	Hydrophilic cream, non-ionic...	7B020413	2
apomix GmbH/Phönix	Hydrophilic cream, non-ionic...	7B020413	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Hydrophilic cream, non-ionic...	7B020413	1
Pharmachem	Hydrophilic cream, non-ionic...	80611	1
Apomix	Hydrophilic cream, non-ionic...	9B300413	3
apomix	Hydrophilic cream, non-ionic...	9A060213	1
AMH Niemann/K	Hydrophilic cream, non-ionic...	10A190312	1
Apomix	Hydrophilic cream, non-ionic...	10B220513	1
Gehe	Hydrophilic cream, non-ionic...	114761	1
Anzag	Hydrophilic cream, non-ionic...	114761	1
Apomix	Hydrophilic cream, non-ionic...	11A040313	1
Apomix / Phönix	Hydrophilic cream, non-ionic...	11A040313	1
Apomix/Gehe GUC	Hydrophilic cream, non-ionic...	11B220513	1
Fagron	Hydrophilic cream, non-ionic...	12/26-M06	1
Gehe	Hydrophilic cream, non-ionic...	120211	1
Gehe	Hydrophilic cream, non-ionic...	122131	1
Phönix Gotha	Hydrophilic cream, non-ionic...	125011	3
Phönix Gotha	Hydrophilic cream, non-ionic...	123711	1
Leyh Pharma	Hydrophilic cream, non-ionic...	125011	1
Gehe	Hydrophilic cream, non-ionic...	125011	1
Apomix	Hydrophilic cream, non-ionic...	12A040313	3
All.Healthc./Apomix	Hydrophilic cream, non-ionic...	12A040313	1
Apomix/Alliance Healthc...	Hydrophilic cream, non-ionic...	12A040313	1
AMH Niemann GmbH& Co. K...	Hydrophilic cream, non-ionic...	12B030613	2
Anzag, 12.09.2012	Hydrophilic cream, non-ionic...	12G19-M07	1
Fagron	Hydrophilic cream, non-ionic...	12J04-M06	2
Fagron	Hydrophilic cream, non-ionic...	12I26-M05	1
Großhandel	Hydrophilic cream, non-ionic...	12I26-M06	1
Phönix Gotha	Hydrophilic cream, non-ionic...	13065	1
Phönix Gotha	Hydrophilic cream, non-ionic...	130651	3
Gehe	Hydrophilic cream, non-ionic...	131451	1
Phönix Gotha	Hydrophilic cream, non-ionic...	13227701	1
Phönix Gotha	Hydrophilic cream, non-ionic...	132331	2
Phönix	Hydrophilic cream, non-ionic...	132331	1
Gehe	Hydrophilic cream, non-ionic...	132331	1
Gehe	Hydrophilic cream, non-ionic...	133121	1
Leyh Pharma GmbH/Anzag	Hydrophilic cream, non-ionic...	133121	1
Phönix Gotha	Hydrophilic cream, non-ionic...	133911	2
Fagron	Hydrophilic cream, non-ionic...	13A30-M01	1
Fagron	Hydrophilic cream, non-ionic...	13B06-M08	1
Apomix	Hydrophilic cream, non-ionic...	13B180613	1
Fagron	Hydrophilic cream, non-ionic...	13D05-M07	3
Fagron	Hydrophilic cream, non-ionic...	13D05-M12	1
Fagron	Hydrophilic cream, non-ionic...	13D18-M06	2
Fagron	Hydrophilic cream, non-ionic...	13D05-M09	1
Fagron	Hydrophilic cream, non-ionic...	13D05M08	1
Fagron	Hydrophilic cream, non-ionic...	13D18-M09	2
Großhandel	Hydrophilic cream, non-ionic...	13D18-M06	1
Fagron	Hydrophilic cream, non-ionic...	13D18-M10	3
Fagron	Hydrophilic cream, non-ionic...	13K15-M07	1
Caelo	Hydrophilic cream, non-ionic...	13k15-m09	1
Anzag/AMH/16.05.2012	Hydrophilic cream, non-ionic...	14B020212	1
Fagron	Hydrophilic cream, non-ionic...	14D28-M04	1
Phönix/ Apomix	Hydrophilic cream, non-ionic...	14B300511	1
apomix/Phönix	Hydrophilic cream, non-ionic...	15B090713	2
AMH Niemann/G	Hydrophilic cream, non-ionic...	15A140512	1
Apomix	Hydrophilic cream, non-ionic...	16B290713	1
Pharmachem/Phönix	Hydrophilic cream, non-ionic...	170412	1
Fagron	Hydrophilic cream, non-ionic...	19A280313	1
Phönix Gotha	Hydrophilic cream, non-ionic...	20B230913	2
AMH Niemann/G	Hydrophilic cream, non-ionic...	20A090413	1
Apomix	Hydrophilic cream, non-ionic...	18B190612	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Hydrophilic cream, non-ionic...	19B100913	1
Apomix	Hydrophilic cream, non-ionic...	19B260612	1
Phönix Gotha	Hydrophilic cream, non-ionic...	20B260612	1
Apomix	Hydrophilic cream, non-ionic...	21B141013	1
Apomix/Noweda	Hydrophilic cream, non-ionic...	21B141013	1
apomix	Hydrophilic cream, non-ionic...	22B060812	1
Fagron	Hydrophilic cream, non-ionic...	24A110413	1
Henry Lamotte	Hydrophilic cream, non-ionic...	22B281013	2
Fagron	Hydrophilic cream, non-ionic...	24A110712	1
Apomix	Hydrophilic cream, non-ionic...	238141113	1
Apomix	Hydrophilic cream, non-ionic...	24B180912	1
Euro OTC	Hydrophilic cream, non-ionic...	25A120712	3
Apomix/Phönix	Hydrophilic cream, non-ionic...	25B161213	1
Caelo	Hydrophilic cream, non-ionic...	25B240912	1
Apomix	Hydrophilic cream, non-ionic...	26B151012	1
Bombastus	Hydrophilic cream, non-ionic...	271397	1
Gehe	Hydrophilic cream, non-ionic...	271397	1
Bombastus	Hydrophilic cream, non-ionic...	273128	1
Bombastus	Hydrophilic cream, non-ionic...	273126	1
Gehe	Hydrophilic cream, non-ionic...	274101	1
Bombastus	Hydrophilic cream, non-ionic...	276595	1
Bombastus	Hydrophilic cream, non-ionic...	276598	2
Bombastus	Hydrophilic cream, non-ionic...	276597	2
Bombastus	Hydrophilic cream, non-ionic...	276599	2
Gehe	Hydrophilic cream, non-ionic...	276599	1
Bombastus	Hydrophilic cream, non-ionic...	276600	1
Bombastus	Hydrophilic cream, non-ionic...	276725	2
Bombastus	Hydrophilic cream, non-ionic...	277871	1
Gehe	Hydrophilic cream, non-ionic...	277872	1
Bombastus	Hydrophilic cream, non-ionic...	277877	2
Caelo	Hydrophilic cream, non-ionic...	277877	1
Caelo	Hydrophilic cream, non-ionic...	277878	1
Bombastus	Hydrophilic cream, non-ionic...	277880	1
Caelo	Hydrophilic cream, non-ionic...	277880	1
Bombastus	Hydrophilic cream, non-ionic...	278025	5
Bombastus	Hydrophilic cream, non-ionic...	277881	1
Bombastus	Hydrophilic cream, non-ionic...	278026	2
Bombastus	Hydrophilic cream, non-ionic...	278028	5
Bombastus	Hydrophilic cream, non-ionic...	278029	1
Bombastus	Hydrophilic cream, non-ionic...	279656	2
Bombastus	Hydrophilic cream, non-ionic...	279991	1
AMH Niemann GmbH & Co. KG	Hydrophilic cream, non-ionic...	27A150513	1
apomix	Hydrophilic cream, non-ionic...	27A150513	1
Bombastus	Hydrophilic cream, non-ionic...	280986	2
Bombastus	Hydrophilic cream, non-ionic...	280984	1
Bombastus	Hydrophilic cream, non-ionic...	280987	1
Gehe	Hydrophilic cream, non-ionic...	280988	1
Caelo	Hydrophilic cream, non-ionic...	280991	1
Gehe	Hydrophilic cream, non-ionic...	280989	2
Gehe	Hydrophilic cream, non-ionic...	280992	1
Bombastus	Hydrophilic cream, non-ionic...	280993	3
Pharmachem / Anzag	Hydrophilic cream, non-ionic...	282212	2
Bombastus	Hydrophilic cream, non-ionic...	281774	1
Gehe	Hydrophilic cream, non-ionic...	282212	1
Bombastus	Hydrophilic cream, non-ionic...	282213	4
Bombastus	Hydrophilic cream, non-ionic...	287716	1
Apomix	Hydrophilic cream, non-ionic...	28B221012	1
Apomix	Hydrophilic cream, non-ionic...	30A030613	1
apomix	Hydrophilic cream, non-ionic...	30A170912	1
Apomix	Hydrophilic cream, non-ionic...	31A170912	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Hydrophilic cream, non-ionic...	32A090713	2
Apomix	Hydrophilic cream, non-ionic...	33A051211	1
Apomix Anzag	Hydrophilic cream, non-ionic...	33A240713	1
PKH Halle/ Alliance H.	Hydrophilic cream, non-ionic...	33a240713	1
PKH Halle/ Alliance H.	Hydrophilic cream, non-ionic...	33A240713	2
Alliance Healthcare/Apo...	Hydrophilic cream, non-ionic...	33A240713	1
Apomix / Noweda	Hydrophilic cream, non-ionic...	33A240713	1
Apomix	Hydrophilic cream, non-ionic...	33B271212	1
AMH Niemann GmbH & Co. KG	Hydrophilic cream, non-ionic...	35A070813	1
Apomix / Phönix	Hydrophilic cream, non-ionic...	35A070813	1
Apomix	Hydrophilic cream, non-ionic...	38A051112	1
apomix/GEHE/02.07.2013	Hydrophilic cream, non-ionic...	40A261112	1
Apomix	Hydrophilic cream, non-ionic...	40A281013	1
Fagron	Hydrophilic cream, non-ionic...	412062-3-P	1
Apomix	Hydrophilic cream, non-ionic...	42A191212	1
Apomix	Hydrophilic cream, non-ionic...	54A021213	1
Alliance Healthcare/Apo...	Hydrophilic cream, non-ionic...	55A181213	1
Fagron	Hydrophilic cream, non-ionic...	55A181213	1
Bombastus	Hydrophilic cream, non-ionic...	82216	1
Bombastus	Hydrophilic cream, non-ionic...	96914	1
Caesar &Loretz GmbH	Hydrophilic cream, non-ionic...		3
Eigene	Hydrophilic cream, non-ionic...	R12/0313	1
Fagron	Transdermal base cream	12H17-T02-005161	2
Fagron	Transdermal base cream	1504065	6
Fagron	Transdermal base cream	1507997	4
Fagron	Transdermal base cream	15A28-T10-018922	3
Fagron	Transdermal base cream	15D28-T06-02194	1
Fagron	Transdermal base cream	1607076	3
Caelo	Transdermal base cream	1609813	1
Fagron	Transdermal base cream	1709950/2	1
Fagron	Transdermal base cream	1802648/2	1
Fagron	Transdermal base cream	4315I-07039	1
Fagron	Transdermal base cream	247191-0003	1
Fagron	Transdermal base cream	12F19-U15-273303	1
Fagron	Transdermal base cream	12H17-T02-000281	1
Fagron	Transdermal base cream	13G17-T03-005497	1
Fagron	Transdermal base cream	14B27-T01-010612	3
Fagron	Transdermal base cream	13G17-T03-010523	1
Fagron	Transdermal base cream	14C14-T03-011060	1
Fagro	Transdermal base cream	14J08-T12-016301	2
Fagron	Transdermal base cream	14H19-T07-014889	1
Fagron	Transdermal base cream	14K21-T03-017364	1
Fagron	Transdermal base cream	15A28-T08-018919	1
Fagron	Transdermal base cream	14C14-T03	1
Fagron	Transdermal base cream	14B27-T010612	1
Fagron	Transdermal base cream	14B27-To1-010612	1
Fagron	Transdermal base cream	15D28-T06-021940	1
Fagron	Transdermal base cream	1602049	1
Apomix	Hydrophilic cream, non-ionic...	1A080113	2
Apomix/Kehr	Hydrophilic cream, non-ionic...	2A090215	1
Apomix	Anionic hydrophlic cream wit...	8B100614	1
Apomix/Kehr	Anionic hydrophilic cream SR...	5A050215	1
Caelo	Anionic hydrophilic cream DAB	13168201	2
Caelo	Anionic hydrophilic cream DAB	1315Q-02540	1
Caelo	Anionic hydrophilic cream DAB	13314101	9
Fagron	Base cream, hydrophobic DAC	13B04-t17-001037	1
Caesar & Lorenz	Base cream, hydrophobic DAC	14028002	1
Caesar & Lorenz	Base cream, hydrophobic DAC	14061501	1
AHD 12.5.14, Preis 6,00	Anionic hydrophilic cream DAB		1
VDL;2.7.14	Anionic hydrophilic cream DAB	14118401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	14138701	3
Caelo	Anionic hydrophilic cream DAB	14186101	5
Caelo	Anionic hydrophilic cream DAB	14186501	15
Caelo	Anionic hydrophilic cream DAB	14186801	7
Caelo	Anionic hydrophilic cream wi...	14214001	7
Caelo	Anionic hydrophilic cream wi...	14214301	2
Caelo	Anionic hydrophilic cream wi...	14214201	32
pharmachem/spangro	Base cream, hydrophobic DAC	142154	1
Caelo	Anionic hydrophilic cream DAB	14233301	5
Caelo	Anionic hydrophilic cream DAB	14285801	5
Caelo	Anionic hydrophilic cream wi...	14334201	12
Caelo	Anionic hydrophilic cream DAB	14341801	1
Caelo	Anionic hydrophilic cream DAB	14341601	7
Caelo	Anionic hydrophilic cream DAB	14374801	3
Caelo	Anionic hydrophilic cream DAB	15001201	9
Sanacorp	Hydrophilic cream, non-ionic...	15039001	1
Caelo	Anionic hydrophilic cream wi...	15048401	6
Caelo	Anionic hydrophilic cream wi...	15048501	16
Caelo	Anionic hydrophilic cream wi...	15048502	12
Caelo	Anionic hydrophilic cream wi...	15132701	5
Caelo	Anionic hydrophilic cream wi...	15132902	14
Caelo	Anionic hydrophilic cream wi...	15133001	13
Caelo	Anionic hydrophilic cream DAB	15200001	4
Caelo	Anionic hydrophilic cream DAB	15267401	14
Caelo	Anionic hydrophilic cream wi...	15273001	28
Caelo	Anionic hydrophilic cream wi...	15273301	11
Caelo	Anionic hydrophilic cream DAB	15400601	1
Caelo	Anionic hydrophilic cream DAB	15397101	3
Caelo	Anionic hydrophilic cream DAB	15401601	1
Caelo	Anionic hydrophilic cream wi...	15405401	18
Caelo	Anionic hydrophilic cream wi...	15405201	5
Fagron	Anionic hydrophilic cream DAB	16046701	1
Caelo	Anionic hydrophilic cream DAB	16046701	1
sanacorp, 5,75EUR, 12.09....	Anionic hydrophilic cream wi...	16191601	2
Caelo	Anionic hydrophilic cream DAB	162500001	1
Sanacorp, 16.03.17, 6,00EUR	Anionic hydrophilic cream DAB	16264501	1
Bombastus	Anionic hydrophilic cream SR...	17001700	1
Bombastus	Anionic hydrophilic cream SR...	18001476	1
Pharmachem/AHD	Base cream, hydrophobic DAC	231213	2
Bombastus	Anionic hydrophilic cream DAB	283480	3
Bombastus	Anionic hydrophilic cream DAB	283483	1
Bombastus	Anionic hydrophilic cream DAB	285589	1
Bombastus	Anionic hydrophilic cream DAB	285593	1
Bombastus	Hydrophilic cream, non-ionic...	286325	1
Bombastus	Anionic hydrophilic cream DAB	287435	1
Bombastus	Anionic hydrophilic cream DAB	287560	1
Bombastus	Base cream, hydrophobic DAC	289717	1
Apomix/Phönix	Anionic hydrophilic cream SR...	28A040815	1
Bombastus	Hydrophilic cream, non-ionic...	291387	1
Bombastus	Anionic hydrophilic cream DAB	290149	2
Bombastus	Anionic hydrophilic cream DAB	291655	1
Bombastus	Anionic hydrophilic cream DAB	291832	2
Sanacorp;21.04.15;9,31EUR	Anionic hydrophilic cream SR...	292977	1
AHD	Anionic hydrophilic cream SR...	292976	1
AHD	Anionic hydrophilic cream SR...	292977	1
Bombastus	Anionic hydrophilic cream DAB	293358	1
Bombastus	Anionic hydrophilic cream SR...	296410	3
AHD	Anionic hydrophilic cream SR...	296409	1
Gehe, 9,99EUR, 29.08.2016	Anionic hydrophilic cream SR...	299255	1
Bombastus	Anionic hydrophilic cream SR...	29926	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream SR...	300954	1
Fagron	Hydrophilic cream, non-ionic...	3713U-03195	1
Fagron	Anionic hydrophilic cream SR...	32A150915	1
AMH Niemann / apomix	Hydrophilic cream, non-ionic...	43a140714	1
Bombastus	Anionic hydrophilic cream SR...	308219	1
Apomix/Phönix	Anionic hydrophilic cream SR...	38A171115	1
Bombastus	Anionic hydrophilic cream SR...	3213I-03327	1
Pharmachem/AHD	Base cream, hydrophobic DAC	110413	2
Sanacorp	Base cream, hydrophobic DAC	13085701	2
Sanacorp	Base cream, hydrophobic DAC	13175902	2
Sanacorp	Base cream, hydrophobic DAC	13275301	2
Euro OTC	Base cream, hydrophobic DAC		1
Pharmachem/AHD	Base cream, hydrophobic DAC	230713	1
Pharmachem/Phoenix	Base cream, hydrophobic DAC	240113	1
Pharmachem/AHD	Base cream, hydrophobic DAC	261112	1
Pharmachem/AHD	Base cream, hydrophobic DAC	26112	1
Gehe	Anionic hydrophilic cream DAB	262935	1
Bombastus	Anionic hydrophilic cream SR...	263823	1
Bombastus	Hydrophilic cream, non-ionic...	265335	1
Bombastus	Hydrophilic cream, non-ionic...	268354	1
Bombastus	Anionic hydrophilic cream SR...	269230	1
Bombastus	Anionic hydrophilic cream DAB	269441	1
Fiebig	Anionic hydrophilic cream DAB	269445	1
Caelo	Anionic hydrophilic cream DAB	270943	2
Bombastus	Anionic hydrophilic cream DAB	271014	2
Gehe	Anionic hydrophilic cream SR...	271701	1
Gehe	Hydrophilic cream, non-ionic...	274216	1
Holdermann	Hydrophilic cream, non-ionic...	1B140113	1
Caelo	Anionic hydrophilic cream DAB	11253201	1
apomix	Anionic hydrophilic cream SR...	5A300113	1
Caelo	Anionic hydrophilic cream DAB	11326601	1
Anzag	Anionic hydrophilic cream DAB	12204601	2
Anzag	Anionic hydrophilic cream DAB	12218501	1
Caelo	Anionic hydrophilic cream DAB	12265001	4
Caelo	Anionic hydrophilic cream DAB	12245801	2
Caelo	Anionic hydrophilic cream DAB	12313501	6
Caelo	Anionic hydrophilic cream DAB	1231401	1
Sanacorp	Anionic hydrophilic cream DAB	12314401	1
Caelo	Anionic hydrophilic cream DAB	12314404	1
Caelo	Anionic hydrophilic cream wi...	13072401	6
Sanacorp	Anionic hydrophilic cream DAB	13154301	1
Caelo	Anionic hydrophilic cream DAB	13154301	6
Caelo	Anionic hydrophilic cream DAB	13167401	2
Caelo	Anionic hydrophilic cream DAB	13167501-3423310192	1
Caelo	Anionic hydrophilic cream DAB	13167801	2
Caelo	Anionic hydrophilic cream DAB	13167501-3439311151	1
Sanacorp	Anionic hydrophilic cream DAB	13168701	1
Caelo	Anionic hydrophilic cream DAB	13169201	5
Caelo	Anionic hydrophilic cream wi...	13244402	10
Caelo	Anionic hydrophilic cream wi...	13244401-3471312131	1
Caelo	Anionic hydrophilic cream DAB	13280201-3452311172	1
Caelo	Anionic hydrophilic cream DAB	13312901	4
Sanacorp	Anionic hydrophilic cream DAB	13313201	1
Caelo	Anionic hydrophilic cream DAB	1331402	1
Caelo	Anionic hydrophilic cream DAB	13364301	4
Caelo	Anionic hydrophilic cream wi...	14030401	10
Caelo	Anionic hydrophilic cream DAB	13314302	4
Caelo	Anionic hydrophilic cream DAB	13409501	2
Caelo	Anionic hydrophilic cream DAB	13346101	4
Anzag	Anionic hydrophilic cream DAB	271710	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream SR...	272513	1
Bombastus	Hydrophilic cream, non-ionic...	274102	1
Caelo	Anionic hydrophilic cream DAB	276319	1
Bombastus	Anionic hydrophilic cream DAB	276316	1
Sanacorp	Anionic hydrophilic cream SR...	274505	1
Gehe	Anionic hydrophilic cream DAB	274765	1
Caelo	Anionic hydrophilic cream DAB	276409	1
Bombastus	Anionic hydrophilic cream DAB	277292	1
Sanacorp	Anionic hydrophilic cream SR...	277566	1
Sanacorp	Anionic hydrophilic cream SR...	277569	1
Bombastus	Anionic hydrophilic cream DAB	278849	1
EUro/AHD	Hydrophilic cream, non-ionic...	277577	1
Gehe	Anionic hydrophilic cream DAB	279057	1
Sanacorp	Anionic hydrophilic cream SR...	279202	1
Bomba/AHD	Anionic hydrophilic cream DAB	280410	1
Bombastus	Anionic hydrophilic cream SR...	279202	1
Bombastus	Anionic hydrophilic cream DAB	280416	1
Sanacorp	Anionic hydrophilic cream SR...	281742	1
Bombastus	Anionic hydrophilic cream SR...	281743	2
Sanacorp	Anionic hydrophilic cream SR...	281809	1
Bombastus	Anionic hydrophilic cream DAB	281998	1
Wolff/AHD	Anionic hydrophilic cream DAB	282107	1
Gehe	Anionic hydrophilic cream DAB	288653	1
Gehe	Hydrophilic cream, non-ionic...	288882	1
Gehe	Anionic hydrophilic cream DAB	290080	1
Apomix	Hydrophilic cream, non-ionic...	30B261112	1
EUro/AHD	Hydrophilic cream, non-ionic...	34A240713	1
AHD	Anionic hydrophilic cream SR...	281866	1
Caelo	Anionic hydrophilic cream DAB	4051506	1
Caelo	Anionic hydrophilic cream DAB	5031508	1
Caelo	Anionic hydrophilic cream DAB	6061501	1
Caelo	Anionic hydrophilic cream DAB	9041505	1
Caelo	Anionic hydrophilic cream DAB	104/AKD2	1
Sanacorp	Anionic hydrophilic cream DAB	104/AKD3	1
Caelo	Anionic hydrophilic cream DAB	12081405	1
Caelo	Anionic hydrophilic cream wi...	12221003	3
Caelo	Anionic hydrophilic cream wi...	13172202	7
Caelo	Anionic hydrophilic cream DAB	13314701	2
Caelo	Anionic hydrophilic cream DAB	133312901	1
Caelo	Anionic hydrophilic cream wi...	13343301	15
Caelo	Anionic hydrophilic cream wi...	13343401	22
Noweda	Anionic hydrophilic cream wi...	13343402	1
Caelo	Anionic hydrophilic cream wi...	13343402	7
Caelo	Anionic hydrophilic cream wi...	13343501	6
Sanacorp	Anionic hydrophilic cream wi...	13343501	3
Caelo	Anionic hydrophilic cream DAB	13346501	4
Caelo	Anionic hydrophilic cream DAB	13363301	6
Alliance Healthcare	Anionic hydrophilic cream DAB	13363401	1
Caelo	Anionic hydrophilic cream DAB	13363401	2
Sanacorp	Anionic hydrophilic cream DAB	13364001	2
Sanacorp 25.02.2014 Ek ...	Anionic hydrophilic cream DAB	13364201	2
Caelo	Anionic hydrophilic cream DAB	13364201	5
Caelo/Sanacorp	Anionic hydrophilic cream DAB	13364201	1
Gehe am 21.02.14	Anionic hydrophilic cream DAB	13364401	1
Caelo	Anionic hydrophilic cream DAB	14016602	2
Caelo	Anionic hydrophilic cream DAB	14016601	1
Fagron	Anionic hydrophilic cream DAB	14016602	1
Caelo	Anionic hydrophilic cream DAB	14016801	5
Caelo	Anionic hydrophilic cream DAB	14017001	2
Caelo	Anionic hydrophilic cream wi...	14030501	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix	Anionic hydrophilic cream wi...	14030601	5
Caelo	Anionic hydrophilic cream wi...	14030601	1
Caelo	Anionic hydrophilic cream wi...	14030701	29
Euro OTC	Anionic hydrophilic cream wi...	14030701	1
Bombastus	Anionic hydrophilic cream wi...	14030701	1
Caelo	Anionic hydrophilic cream wi...	14030702	5
Phoenix	Anionic hydrophilic cream wi...	14030702	1
Caelo	Anionic hydrophilic cream DAB	14042301	1
Großhandel	Anionic hydrophilic cream DAB	14042401	1
Caelo	Anionic hydrophilic cream DAB	14042401	1
Caelo	Anionic hydrophilic cream DAB	14042501	1
Caelo	Anionic hydrophilic cream DAB	14042601	2
PHOENIX 22.05.2014	Anionic hydrophilic cream DAB	14061801	1
GEHE 22.03.2014	Anionic hydrophilic cream DAB	14062301	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	14062401	1
Caelo	Anionic hydrophilic cream DAB	14062501	6
Gehe	Anionic hydrophilic cream DAB	14062501	1
Gehe 01.08.14	Anionic hydrophilic cream DAB	14062501	1
Noweda	Anionic hydrophilic cream DAB	14062701	1
Caelo	Anionic hydrophilic cream DAB	14070801	1
Sanacorp 11.07.14 EK 5,76	Anionic hydrophilic cream DAB	141118401	1
Caelo	Anionic hydrophilic cream wi...	141146101	1
Phoenix	Anionic hydrophilic cream DAB	14118201	1
Bombastus	Anionic hydrophilic cream DAB	14118301	3
ahd/lamotte	Anionic hydrophilic cream DAB	14118401	1
Phoenix 11.07.2014	Anionic hydrophilic cream DAB	14139001	2
Caesar & Loretz	Anionic hydrophilic cream DAB	14144301	1
AHDE 12.09.2014	Anionic hydrophilic cream DAB	14144301	1
Caelo	Anionic hydrophilic cream DAB	14144401	6
Caelo	Anionic hydrophilic cream DAB	14144601	1
Caelo	Anionic hydrophilic cream wi...	14146002	14
Phoenix	Anionic hydrophilic cream wi...	14146001	2
Caelo	Anionic hydrophilic cream wi...	14146003	7
Caelo	Anionic hydrophilic cream DAB	141460201	1
Caelo	Anionic hydrophilic cream wi...	14146101	12
Caelo	Anionic hydrophilic cream DAB	14146501	2
Caelo	Anionic hydrophilic cream DAB	14146601	3
Caelo	Anionic hydrophilic cream DAB	1415034401	1
Caelo	Anionic hydrophilic cream wi...	1416101	1
Caelo	Anionic hydrophilic cream DAB	1418101	2
Caelo	Anionic hydrophilic cream DAB	14185201	5
Caelo	Anionic hydrophilic cream DAB	14185301	10
Caelo	Anionic hydrophilic cream DAB	14185601	7
Caelo	Anionic hydrophilic cream DAB	14185702	1
Caelo	Anionic hydrophilic cream DAB	14185901	4
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	14185801	1
Caelo	Anionic hydrophilic cream DAB	14186001	2
Caelo	Anionic hydrophilic cream DAB	1418601	1
Sanacorp	Anionic hydrophilic cream DAB	14186702	1
Sanacorp	Anionic hydrophilic cream DAB	142/AKD52	1
Caelo	Anionic hydrophilic cream wi...	14214101	8
aelo	Anionic hydrophilic cream wi...	14214101	1
Caelo	Anionic hydrophilic cream wi...	14214201	1
Caelo	Anionic hydrophilic cream wi...	14214202	6
Sanacorp	Anionic hydrophilic cream wi...	14214301	7
Caelo	Anionic hydrophilic cream DAB	14232801	2
Gehe	Anionic hydrophilic cream DAB	14232901	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	14232801	1
Caelo	Anionic hydrophilic cream DAB	14232901	4
Alliance health	Anionic hydrophilic cream DAB	14233001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	1423301	2
Caelo	Anionic hydrophilic cream DAB	142331001	1
sanacorp	Anionic hydrophilic cream DAB	14233101	1
Sana 20.11.14 EK 5,76	Anionic hydrophilic cream DAB	14233101	1
Alliance21.01.15	Anionic hydrophilic cream DAB	14233101	1
Phönix	Anionic hydrophilic cream DAB	14233401	1
Caelo	Anionic hydrophilic cream DAB	14233501	2
Caelo	Anionic hydrophilic cream DAB	14255401	1
Caelo	Anionic hydrophilic cream wi...	1424202	1
AHD 05.03.15 3,33	Anionic hydrophilic cream DAB	14264001	1
DAB 2008	Anionic hydrophilic cream DAB	14264001	1
Hedinger/Noweda	Anionic hydrophilic cream DAB	14264001	1
Phönix	Anionic hydrophilic cream DAB	14264301	1
Caelo	Anionic hydrophilic cream DAB	14264302	1
Sanacorp 19.12.14 EK 5,76	Anionic hydrophilic cream DAB	14285601	1
Sanacorp	Anionic hydrophilic cream DAB	14285601	1
Caelo	Anionic hydrophilic cream DAB	14285701	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	14294201	1
Caelo	Anionic hydrophilic cream DAB	14294301	2
Caelo	Anionic hydrophilic cream DAB	14294601	2
Caelo	Anionic hydrophilic cream DAB	14294701	3
Caelo	Anionic hydrophilic cream DAB	14294801	2
Phoenix 14.01.2015 Ek: ...	Anionic hydrophilic cream DAB	14295101	1
Caelo	Anionic hydrophilic cream DAB	14295101	3
Caelo	Anionic hydrophilic cream DAB	14295201	3
Caelo	Anionic hydrophilic cream wi...	14334301	24
Sanacorp	Anionic hydrophilic cream wi...	14334401	1
Euro OTC	Anionic hydrophilic cream wi...	14334501	1
Caelo	Anionic hydrophilic cream wi...	14334502	6
sanacorp	Anionic hydrophilic cream DAB	14341601	1
Caelo	Anionic hydrophilic cream DAB	14341701	8
Caelo	Anionic hydrophilic cream DAB	14342001	1
Caelo	Anionic hydrophilic cream DAB	14342002	2
Caelo	Anionic hydrophilic cream DAB	14342101	9
Caelo	Anionic hydrophilic cream DAB	14342901	2
WE:16.04.2016 Sanacor...	Anionic hydrophilic cream DAB	14342901	1
Caelo	Anionic hydrophilic cream DAB	14366802	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	15027301	1
caesar&loretz,Phönix	Anionic hydrophilic cream DAB	15034301	1
Sanacorp	Anionic hydrophilic cream DAB	15034401	1
Caelo	Anionic hydrophilic cream DAB	15034501	2
Caelo	Anionic hydrophilic cream wi...	15048301	15
Caelo	Anionic hydrophilic cream wi...	1512901	1
Caelo	Anionic hydrophilic cream wi...	15132801	9
Caelo	Anionic hydrophilic cream wi...	15132901	4
Caelo	Anionic hydrophilic cream wi...	1513301	3
Caelo	Anionic hydrophilic cream DAB	15199901	6
Sanacorp	Anionic hydrophilic cream DAB	15199901	1
lamotte/gehe	Anionic hydrophilic cream DAB	15200001	1
Caelo	Anionic hydrophilic cream DAB	15200202	1
Caelo	Anionic hydrophilic cream DAB	15200401	1
Cealo/Anzag	Anionic hydrophilic cream DAB	15200401	1
Caesar&loretz,Anzag	Anionic hydrophilic cream DAB	15258901	1
Caelo	Anionic hydrophilic cream DAB	15259301	1
Caelo	Anionic hydrophilic cream DAB	15267001	2
phönix	Anionic hydrophilic cream DAB	15267401	1
Caelo	Anionic hydrophilic cream DAB	15267701	3
Caelo	Anionic hydrophilic cream wi...	15272901	2
Caelo	Anionic hydrophilic cream wi...	15273101	13
Caelo	Anionic hydrophilic cream wi...	1527301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Anionic hydrophilic cream wi...	15273101	1
Caelo	Anionic hydrophilic cream wi...	15273201	14
Rondell Apotheke	Anionic hydrophilic cream wi...	15273301	1
Fischer/Gehe	Anionic hydrophilic cream DAB	15315601	1
Caelo	Anionic hydrophilic cream DAB	1535101	1
Caelo	Anionic hydrophilic cream DAB	153552901	1
Caelo	Anionic hydrophilic cream DAB	15357101	1
Caelo	Anionic hydrophilic cream DAB	15357201	3
Caelo	Anionic hydrophilic cream DAB	15357501	2
Caelo	Anionic hydrophilic cream DAB	15396901	1
Caelo	Anionic hydrophilic cream DAB	15400501	2
Noweda	Anionic hydrophilic cream DAB	15400501	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	15400601	1
Sanacorp 21,01.16	Anionic hydrophilic cream DAB	15400901	1
Sanacorp	Anionic hydrophilic cream DAB	15401201	1
Caelo	Anionic hydrophilic cream DAB	15401301	2
Caelo	Anionic hydrophilic cream DAB	15401302	4
Caelo	Anionic hydrophilic cream DAB	15401501	2
Caelo	Anionic hydrophilic cream wi...	1540521	1
Caeol/AHD	Anionic hydrophilic cream wi...	15405301	1
Caelo	Anionic hydrophilic cream wi...	15405301	8
Caelo	Anionic hydrophilic cream DAB	1542801	1
Caelo	Anionic hydrophilic cream DAB	15428301	7
Sana 23.03.16 Ek 5,73	Anionic hydrophilic cream DAB	15428401	1
Caelo	Anionic hydrophilic cream DAB	15431901	2
Caelo	Anionic hydrophilic cream DAB	15431801	3
Caelo	Anionic hydrophilic cream wi...	1548301	2
Caelo	Anionic hydrophilic cream DAB	16046501	1
Fagron	Anionic hydrophilic cream DAB	16047301	1
Caelo	Anionic hydrophilic cream DAB	16058101	2
Gehe	Anionic hydrophilic cream wi...	16067301	1
Sanacorp	Anionic hydrophilic cream wi...	16067901	1
Caelo	Anionic hydrophilic cream DAB	16071506	1
Caelo	Anionic hydrophilic cream DAB	161431	1
Caelo	Anionic hydrophilic cream DAB	16150601	1
Caelo	Anionic hydrophilic cream DAB	16156201	2
Caelo	Anionic hydrophilic cream DAB	16250102	1
Caelo	Anionic hydrophilic cream DAB	16258901	1
Caelo	Anionic hydrophilic cream DAB	17101501	1
Caelo	Anionic hydrophilic cream wi...	17030310	1
Caelo	Anionic hydrophilic cream wi...	17156602	1
Caelo	Anionic hydrophilic cream DAB	2111401	1
Caelo	Anionic hydrophilic cream wi...	20071504	1
Caelo	Anionic hydrophilic cream wi...	2114M25401	1
Caelo	Anionic hydrophilic cream DAB	2414E2540	1
Sanacorp	Anionic hydrophilic cream wi...	24AKD1	1
Caelo	Anionic hydrophilic cream wi...	2514U25401	1
Caelo	Anionic hydrophilic cream DAB	27041512	1
Caelo	Anionic hydrophilic cream wi...	27031504	1
Caesar & Loretz GmbH Hi...	Anionic hydrophilic cream DAB	280291	1
	Anionic hydrophilic cream DAB	285597	1
sanacorp	Anionic hydrophilic cream DAB	286501	1
Bombastus	Anionic hydrophilic cream DAB	288575	1
sanacorp	Anionic hydrophilic cream DAB	288655	1
Bombastus	Anionic hydrophilic cream DAB	288720	1
Caelo	Anionic hydrophilic cream DAB	290145	1
Bombastus	Anionic hydrophilic cream DAB	290075	1
Bombastus	Anionic hydrophilic cream DAB	290218	2
Bombastus	Anionic hydrophilic cream DAB	293361	1
	Anionic hydrophilic cream DAB	293437	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Anionic hydrophilic cream DAB	293446	1
Bombastus	Anionic hydrophilic cream DAB	295399	2
Bombastus	Anionic hydrophilic cream DAB	295404	1
Bombastus	Anionic hydrophilic cream DAB	295401	1
Bombastus	Anionic hydrophilic cream DAB	296704	1
Bombastus	Anionic hydrophilic cream DAB	296781	1
Bombastus	Anionic hydrophilic cream DAB	297666	1
Bombastus	Anionic hydrophilic cream DAB	297673	1
Caelo	Anionic hydrophilic cream DAB	300353	1
Caelo	Anionic hydrophilic cream DAB	30091403	1
Bombastus	Anionic hydrophilic cream DAB	301266	1
Bombastus	Anionic hydrophilic cream DAB	302465	1
Caelo	Anionic hydrophilic cream DAB	304105	1
Bombastus	Anionic hydrophilic cream DAB	305383	1
Bombastus	Anionic hydrophilic cream DAB	307039	1
Caelo	Anionic hydrophilic cream wi...	32/AKD2	1
Sanacorp	Anionic hydrophilic cream wi...	39/AKD4	1
Sanacorp	Anionic hydrophilic cream wi...	41/AKD6	1
Caelo	Anionic hydrophilic cream wi...	45/AKD7	1
Caelo	Anionic hydrophilic cream wi...	54/AKD12	1
Sanacorp	Anionic hydrophilic cream wi...	50/AKD12	1
Sanacorp	Anionic hydrophilic cream wi...		1
Caelo	Anionic hydrophilic cream wi...	CHB13343401	1
Caesar & Loretz GmbH Hi...	Anionic hydrophilic cream wi...	B1334301	1
Caelo	Anionic hydrophilic cream wi...	112210016	1
Caelo	Anionic hydrophilic cream DAB	112212540	1
Phönix	Anionic hydrophilic cream DAB	13265602	1
eigene Herstellung	Anionic hydrophilic cream DAB	2/270114	1
Caelo	Anionic hydrophilic cream DAB	314Q2540	1
Caelo	Anionic hydrophilic cream DAB	6021405	1
Caelo	Anionic hydrophilic cream DAB	6021406	1
Caelo	Anionic hydrophilic cream wi...	7011403	1
Caelo	Anionic hydrophilic cream DAB	8051301	1
Caelo	Anionic hydrophilic cream DAB	8111301	1
Caelo	Anionic hydrophilic cream DAB	714M2540	1
Caelo	Anionic hydrophilic cream SR...	9101202	1
Caelo	Anionic hydrophilic cream DAB	10092013A	1
eigene Herstellung	Anionic hydrophilic cream DAB	100214-4	1
Caelo	Anionic hydrophilic cream DAB	1108053-02	1
Caelo	Anionic hydrophilic cream DAB	1114U2540	1
Caelo	Anionic hydrophilic cream SR...	1109212540	1
Gehe	Anionic hydrophilic cream DAB	11219901	1
noweda	Anionic hydrophilic cream DAB	11238910	1
Caelo	Anionic hydrophilic cream DAB	113313001	2
Caelo	Anionic hydrophilic cream DAB	12010501	1
Caelo	Anionic hydrophilic cream DAB	12010601	1
Caelo	Anionic hydrophilic cream DAB	12026001	1
Caelo	Anionic hydrophilic cream DAB	12084101	3
Caelo	Anionic hydrophilic cream DAB	12100101	1
Caesar & Loretz GmbH, 4...	Anionic hydrophilic cream DAB	12100601	1
Caelo	Anionic hydrophilic cream DAB	12101101	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	12101101	1
Caelo	Anionic hydrophilic cream DAB	12101302	1
Caelo	Anionic hydrophilic cream wi...	1213M-07009	1
Caelo	Anionic hydrophilic cream DAB	1218501	1
Noweda	Anionic hydrophilic cream DAB	1218701	1
Caelo	Anionic hydrophilic cream DAB	12188001	1
Caelo	Anionic hydrophilic cream DAB	12187801	2
Caelo	Anionic hydrophilic cream DAB	12188101	1
Caelo	Anionic hydrophilic cream DAB	12195301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream wi...	12216501	1
Noweda	Anionic hydrophilic cream DAB	12204701	1
Caelo	Anionic hydrophilic cream DAB	12217801	5
Noweda	Anionic hydrophilic cream DAB	12217801	1
Caelo	Anionic hydrophilic cream DAB	12218001	2
Caelo	Anionic hydrophilic cream DAB	12217901	2
Phönix 12.10.2012	Anionic hydrophilic cream DAB	12218201	1
Caelo	Anionic hydrophilic cream DAB	12218201	2
Sanacorp	Anionic hydrophilic cream DAB	12218501	2
Caelo	Anionic hydrophilic cream DAB	12218701	1
Caelo	Anionic hydrophilic cream DAB	12218801	3
Caelo	Anionic hydrophilic cream wi...	12221002	5
vdl	Anionic hydrophilic cream wi...	12221002	1
Caelo	Anionic hydrophilic cream wi...	1222102	1
Caelo	Anionic hydrophilic cream DAB	12245503	1
Sanacorp	Anionic hydrophilic cream DAB	12264401	1
Caelo	Anionic hydrophilic cream DAB	12264401	1
Caelo	Anionic hydrophilic cream DAB	12264501	5
Hedinger	Anionic hydrophilic cream DAB	12264501	1
Caelo	Anionic hydrophilic cream DAB	12264701	2
Caelo	Anionic hydrophilic cream DAB	12264502	1
Caelo	Anionic hydrophilic cream DAB	12264901	1
Gehe	Anionic hydrophilic cream DAB	12265001	1
Caelo	Anionic hydrophilic cream DAB	12265101	1
Sanacorp	Anionic hydrophilic cream DAB	12265301	1
Fagron	Anionic hydrophilic cream DAB	12265501	1
Phoenix 25.01.2013	Anionic hydrophilic cream DAB	12265901	1
Caelo	Anionic hydrophilic cream DAB	12265901	1
Caelo	Anionic hydrophilic cream wi...	12283001	1
Caelo	Anionic hydrophilic cream wi...	12301301	2
GEHE	Anionic hydrophilic cream wi...	12301302	1
Caelo	Anionic hydrophilic cream wi...	12301302	5
Fagron	Anionic hydrophilic cream DAB	12313202	1
Caelo	Anionic hydrophilic cream DAB	12313301	4
Gehe	Anionic hydrophilic cream DAB	12313401	1
Caelo	Anionic hydrophilic cream DAB	12313401	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	12313501	1
Post-Apotheke	Anionic hydrophilic cream DAB	12313701	1
Phoenix 26.03.13	Anionic hydrophilic cream DAB	12313901	1
Caelo	Anionic hydrophilic cream DAB	12314201	4
Großhandel	Anionic hydrophilic cream DAB	12314201	1
Noweda	Anionic hydrophilic cream DAB	12314301	1
Sanacorp	Anionic hydrophilic cream DAB	12314301	1
Caelo	Anionic hydrophilic cream DAB	12314501	2
Großhandel	Anionic hydrophilic cream DAB	12314501	1
Phönix 10.04.2013	Anionic hydrophilic cream DAB	12314501	1
Phönix	Anionic hydrophilic cream DAB	1231451	1
Großhandel	Anionic hydrophilic cream DAB	12314901	1
Caelo	Anionic hydrophilic cream DAB	1234301	1
Caelo	Anionic hydrophilic cream wi...	12344401	1
Caelo	Anionic hydrophilic cream DAB	1234569	1
Caelo	Anionic hydrophilic cream wi...	12352701	5
Caelo	Anionic hydrophilic cream wi...	12352702	4
Caesar & Loretz GmbH, 4...	Anionic hydrophilic cream wi...	12352702	1
Caelo	Anionic hydrophilic cream wi...	12352703	2
Caelo	Anionic hydrophilic cream DAB	13035902	4
Noweda	Anionic hydrophilic cream DAB	13035902	1
Phönix	Anionic hydrophilic cream DAB	13036001	1
Caelo	Anionic hydrophilic cream DAB	13036302	1
Phoenix	Anionic hydrophilic cream DAB	13038201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Anionic hydrophilic cream DAB	13038201	1
Caelo	Anionic hydrophilic cream DAB	13071304	1
Caelo	Anionic hydrophilic cream wi...	13072301	2
Caelo	Anionic hydrophilic cream wi...	13072302	12
Caesar & Loretz GmbH, 4...	Anionic hydrophilic cream wi...	13072302	3
Noweda	Anionic hydrophilic cream wi...	13072302	1
Caelo	Anionic hydrophilic cream wi...	13072303	3
GEHE	Anionic hydrophilic cream wi...	13072303	1
Caelo	Anionic hydrophilic cream wi...	130723301	1
noweda	Anionic hydrophilic cream wi...	13072401	1
Noweda	Anionic hydrophilic cream wi...	13072401	1
PHÖNIX 06.06.2013	Anionic hydrophilic cream wi...	13072401	1
Caelo	Anionic hydrophilic cream SR...	1310212540	1
Caelo	Anionic hydrophilic cream wi...	13125901	9
Caesar & Loretz GmbH, 4...	Anionic hydrophilic cream wi...	13125901	2
Caelo	Anionic hydrophilic cream wi...	13125902	3
Caelo	Anionic hydrophilic cream wi...	13125903	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	13128701	1
Sanacorp	Anionic hydrophilic cream DAB	13154101	1
Gehe	Anionic hydrophilic cream DAB	13154201	2
Caelo	Anionic hydrophilic cream DAB	13154401	1
Cealo	Anionic hydrophilic cream wi...	13165001	1
Großhandel	Anionic hydrophilic cream DAB	13165101	1
Gehe	Anionic hydrophilic cream DAB	13165701	1
Noweda	Anionic hydrophilic cream DAB	13167401	1
Caelo	Anionic hydrophilic cream DAB	13167701	2
Caelo	Anionic hydrophilic cream DAB	13168001	7
Sanacorp	Anionic hydrophilic cream DAB	13168001	1
Caelo	Anionic hydrophilic cream DAB	13168301	5
Caelo	Anionic hydrophilic cream DAB	13168501	10
Noweda	Anionic hydrophilic cream DAB	13168601	1
Caelo	Anionic hydrophilic cream DAB	13168701	6
Caelo	Anionic hydrophilic cream DAB	13169101	4
Großhandel	Anionic hydrophilic cream DAB	13169201	1
Caelo	Anionic hydrophilic cream wi...	13172201	1
Phönix	Anionic hydrophilic cream wi...	13172203	1
Noweda	Anionic hydrophilic cream wi...	13172203	1
Caelo	Anionic hydrophilic cream wi...	13172203	4
Caelo	Anionic hydrophilic cream wi...	1324402	1
Caelo	Anionic hydrophilic cream wi...	1324401	1
Phönix/gepepharm	Anionic hydrophilic cream wi...	13244201	2
Caelo	Anionic hydrophilic cream wi...	13244201	9
Caesar & Loretz GmbH, 4...	Anionic hydrophilic cream wi...	13244201	1
Caelo	Anionic hydrophilic cream wi...	13244301	4
Caelo	Anionic hydrophilic cream wi...	13244302	3
Caesar & Loretz GmbH, 4...	Anionic hydrophilic cream wi...	13244401	1
Caelo	Anionic hydrophilic cream wi...	13244401	17
Gehe	Anionic hydrophilic cream wi...	13244401	1
Phoenix 26.11.13	Anionic hydrophilic cream DAB	13265501	1
Caelo	Anionic hydrophilic cream DAB	13265701	1
Noweda	Anionic hydrophilic cream DAB	13270801	1
Caelo	Anionic hydrophilic cream DAB	13271001	3
Cealo/Sanacorp	Anionic hydrophilic cream DAB	13275501	1
Caelo	Anionic hydrophilic cream DAB	13275501	1
Caelo	Anionic hydrophilic cream DAB	13275701	1
Caelo	Anionic hydrophilic cream DAB	13280001	3
Noweda	Anionic hydrophilic cream DAB	13280101	1
Gehe	Anionic hydrophilic cream DAB	13280201	1
Sanacorp	Anionic hydrophilic cream DAB	13280201	1
Caelo	Anionic hydrophilic cream DAB	13310001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Großhandel	Anionic hydrophilic cream DAB	13312901	1
Sanacorp	Anionic hydrophilic cream DAB	13313101	1
AHD, 31.10.2013	Anionic hydrophilic cream DAB	13313201	1
Phönix 11.11.2014	Anionic hydrophilic cream DAB	13313301	1
Caesar & Loretz	Anionic hydrophilic cream DAB	13313301	1
Caelo	Anionic hydrophilic cream DAB	13313401	5
Großhandel	Anionic hydrophilic cream DAB	13313601	1
Caelo	Anionic hydrophilic cream DAB	13314601	2
Phoenix	Anionic hydrophilic cream DAB	13314901	1
Caelo	Anionic hydrophilic cream DAB	13314901	1
Caelo	Anionic hydrophilic cream DAB	1331501	1
Caelo	Anionic hydrophilic cream wi...	1334302	1
Noweda	Anionic hydrophilic cream wi...	13343301	1
Caelo	Anionic hydrophilic cream wi...	13344301	1
Caelo	Anionic hydrophilic cream DAB	13346001	4
Phoenix 24.01.2014	Anionic hydrophilic cream DAB	13346101	1
Phönix	Anionic hydrophilic cream DAB	13346301	1
Großhandel	Anionic hydrophilic cream DAB	13346401	1
Caesar & Loretz GmbH Hi...	Anionic hydrophilic cream DAB	13346401	1
Caelo	Anionic hydrophilic cream DAB	13346401	1
Caelo	Anionic hydrophilic cream DAB	13346701	1
Großhandel	Anionic hydrophilic cream DAB	13346701	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	13363601	1
Großhandel	Anionic hydrophilic cream DAB	1369201	1
Phoenix	Anionic hydrophilic cream DAB	13364201	1
Caelo	Anionic hydrophilic cream DAB	1368001	1
Caesar & Loretz GmbH	Anionic hydrophilic cream DAB	13409401	1
Phoenix	Anionic hydrophilic cream DAB	14016602	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream wi...	14030701	1
Caelo	Anionic hydrophilic cream DAB	14042201	1
Gehe	Anionic hydrophilic cream DAB	14042201	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	14118201	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	14139001	1
Phoenix	Anionic hydrophilic cream DAB	14139001	1
Caesar&lorenz,phönix	Anionic hydrophilic cream DAB	14186702	1
Caelo	Anionic hydrophilic cream DAB	142/AKD51	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	14233601	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	14285601	1
Caesar&lorenz,Phoenix	Anionic hydrophilic cream DAB	14294801	1
Caelo	Anionic hydrophilic cream DAB	1601312540	2
Caelo	Anionic hydrophilic cream DAB	19101301	1
Caelo	Anionic hydrophilic cream DAB	20301012-3	1
Caelo	Anionic hydrophilic cream SR...	21081201	1
Caelo	Anionic hydrophilic cream DAB	20021311	1
Caelo	Anionic hydrophilic cream DAB	22091202	1
Caelo	Anionic hydrophilic cream wi...	18/AKD1511	1
Euro OTC	Anionic hydrophilic cream DAB	237563	1
Caelo	Anionic hydrophilic cream DAB	24021401	1
Caelo	Anionic hydrophilic cream DAB	2413Q02540	1
Bombastus	Anionic hydrophilic cream DAB	273538	1
Anzag, 28.11.12	Anionic hydrophilic cream DAB	274696	1
Bombastus	Anionic hydrophilic cream DAB	274696	1
Bombastus	Anionic hydrophilic cream DAB	274698	1
Phönix	Anionic hydrophilic cream DAB	276433	1
Bombastus	Anionic hydrophilic cream DAB	277462	1
Bombastus	Anionic hydrophilic cream DAB	280288	2
Bombastus	Anionic hydrophilic cream DAB	280413	1
Bombastus	Anionic hydrophilic cream DAB	282107	1
Bombastus	Anionic hydrophilic cream DAB	283412	1
Caelo	Anionic hydrophilic cream wi...	2913M25401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Anionic hydrophilic cream DAB	30081301	1
Caelo	Anionic hydrophilic cream wi...	3125902	1
Caelo	Anionic hydrophilic cream DAB	3313U-02540	1
Caelo	Anionic hydrophilic cream DAB	3413M-02540	1
Caelo	Anionic hydrophilic cream DAB	3321	1
Caelo	Anionic hydrophilic cream DAB	3713A02540	1
Caelo	Anionic hydrophilic cream DAB	3960/V07/15	1
Caelo	Anionic hydrophilic cream DAB	3965/V07/15	1
Caelo	Anionic hydrophilic cream DAB	4513E2540	1
Caelo	Anionic hydrophilic cream DAB	4513U2540	1
Caelo	Anionic hydrophilic cream wi...	Ch13125901	1
Caelo	Anionic hydrophilic cream DAB	Ch13128601	1
Caelo	Anionic hydrophilic cream DAB	CH13314001	1
Caelo	Anionic hydrophilic cream DAB	pzn01119720	1
Caelo	Anionic hydrophilic cream DAB	U-030513-1	1
Caelo	Anionic hydrophilic cream DAB	U-310513-1	1
Dr. Wolff	Wolff [®] base half-fat cream	48882171377	1
Dr. Wolff	Wolff [®] base half-fat cream	12081508	1
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	170420WBC/605551	1
Dr. August Wolff GmbH &...	Wolff [®] base half-fat cream		1
Dr. August Wolff GmbH &...	Wolff [®] base half-fat cream	170519WBC/610991	1
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	170616WBC/610991	1
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	170920WBC/610991	1
Dr. August Wolff GmbH &...	Wolff [®] base half-fat cream	171113WBC/701451	1
Dr. August Wolff GmbH &...	Wolff [®] base half-fat cream	180312WB/705760	1
Dr. Wolff	Wolff [®] base half-fat cream	410400	1
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	434360	4
VDL;16.04.15;46,52EUR	Wolff [®] base half-fat cream	434360	2
Dr. Wolff/Sanacorp	Wolff [®] base half-fat cream	434360	1
Dr. August Wolff	Wolff [®] base half-fat cream	434360	1
Wolff/ Gehe	Wolff [®] base half-fat cream	434360	1
Dr. Wolff	Wolff [®] base half-fat cream	436750	3
Wolff/Alliance	Wolff [®] base half-fat cream	436750	1
Dr. Wolff	Wolff [®] base half-fat cream	438160	1
Wolff/Phönix	Wolff [®] base half-fat cream	438160	1
Dr. Wolff GmbH & Co.KG A...	Wolff [®] base half-fat cream	438160	1
Dr. Wolff / Kehr	Wolff [®] base half-fat cream	438160	1
A. Wolff/Kehr	Wolff [®] base half-fat cream	438340	2
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	438340	1
wolff GmbH/AHD 21.09.2015	Wolff [®] base half-fat cream	438340	1
Wolff/Phoenix	Wolff [®] base half-fat cream	4715E-03719	1
Wolff	Wolff [®] base half-fat cream	501080	3
Wolff GmbH	Wolff [®] base half-fat cream	501080	1
Dr. Wolff	Wolff [®] base half-fat cream	501080	1
Dr. Wolff GmbH	Wolff [®] base half-fat cream	501080	1
Dr. Wolff	Wolff [®] base half-fat cream	501080	1
Dr. Wolff	Wolff [®] base half-fat cream	501080	1
Sanacorp/Dr. Wolff	Wolff [®] base half-fat cream	501082	1
Dr. Wolff	Wolff [®] base half-fat cream	501082	1
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	501082	1
Dr. Wolff	Wolff [®] base half-fat cream	501470	4
Wolff/Phönix	Wolff [®] base half-fat cream	501470	1
Caelo	Wolff [®] base half-fat cream	501470	6
Dr. Wolff/ Phoenix	Wolff [®] base half-fat cream	501470	2
Wolff GmbH / Noweda	Wolff [®] base half-fat cream	501470	1
Wolff/ Anzag	Wolff [®] base half-fat cream	501470	1
Wolff/Sanacorp	Wolff [®] base half-fat cream	501470	1
Dr. Wolff /Gehe	Wolff [®] base half-fat cream	501870	1
Wolff	Wolff [®] base half-fat cream	506360	1
Dr. Wolff/Gehe	Wolff [®] base half-fat cream	506360	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Dr. August Wolff GmbH	Wolff [®] base half-fat cream	506360	1
Chiesi / Sanacorp/EK: 6...	Wolff [®] base half-fat cream	506360	1
Dr. Wolff/Phoenix	Wolff [®] base half-fat cream	506360	1
Caelo	Wolff [®] base half-fat cream	506361	1
Dr. August Wolff GmbH/S...	Wolff [®] base half-fat cream	506361	1
Dr. Wolff	Wolff [®] base half-fat cream	510400	1
Dr. Wollff	Wolff [®] base half-fat cream	510400	1
Euro OTC	Wolff [®] base half-fat cream	510400	1
Gehe/ Wolff	Wolff [®] base half-fat cream	510400	1
Dr Wolff	Wolff [®] base half-fat cream	510400	1
Wolff/AHD	Wolff [®] base half-fat cream	510621	1
Dr. Wolff/Phönix	Wolff [®] base half-fat cream	510621	1
Dr. Wolff / Kehr	Wolff [®] base half-fat cream	510621	1
Wolf	Wolff [®] base half-fat cream	511420	1
Dr Wolff/Noweda	Wolff [®] base half-fat cream	511420	3
Caelo	Wolff [®] base half-fat cream	511420	2
Dr.A. Wolff, Sanacorp	Wolff [®] base half-fat cream	511420	2
Dr. August Wolff GmbH	Wolff [®] base half-fat cream	511420	2
Dr. August Wolff	Wolff [®] base half-fat cream	511420	1
Dr. Wolff, Sanacorp ,45,7...	Wolff [®] base half-fat cream	511420	3
Phoenix	Wolff [®] base half-fat cream	511420	1
Wolff/Alliance	Wolff [®] base half-fat cream	511420	1
Wolff/Noweda	Wolff [®] base half-fat cream	511420	1
WE: Sanacorp 18.10.2016...	Wolff [®] base half-fat cream	601880	1
Chiesi	Wolff [®] base half-fat cream	601880	1
Dr.A. Wolff	Wolff [®] base half-fat cream	601880	1
Wolff/Sanacorp	Wolff [®] base half-fat cream	602050	3
Caelo	Wolff [®] base half-fat cream	602050	1
Caelo	Wolff [®] base half-fat cream	601880	1
Phoenix	Wolff [®] base half-fat cream	602050	3
sanacorp, 43,56EUR , 10.0...	Wolff [®] base half-fat cream	602050	4
Dr. August Wolff GmbH	Wolff [®] base half-fat cream	602050	2
Taoasis	Wolff [®] base half-fat cream	602050	1
Dr. August Wolff	Wolff [®] base half-fat cream	602050	3
Caelo	Wolff [®] base half-fat cream	602290	1
AHD	Wolff [®] base half-fat cream	602290	1
Dr.A. Wolff/Gehe	Wolff [®] base half-fat cream	602290	1
Dr. Wolff / Kehr	Wolff [®] base half-fat cream	602290	1
Dr. Wolff	Wolff [®] base half-fat cream	605551	1
Dr. Wolff / Kehr	Wolff [®] base half-fat cream	605700	3
Noweda/Wolff GmbH	Wolff [®] base half-fat cream	602290	1
Dr. August Wolff	Wolff [®] base half-fat cream	605700	1
Dr. Wolff/AHD	Wolff [®] base half-fat cream	606050	1
Dr Wolff	Wolff [®] base half-fat cream	606050	3
Caelo	Wolff [®] base half-fat cream	606050	1
Sanacorp, 11.03.17, 45,74EUR	Wolff [®] base half-fat cream	608180	2
Dr. Wolff/Sanacorp	Wolff [®] base half-fat cream	608180	3
Dr. Wolff / OG	Wolff [®] base half-fat cream	608180	1
Dr. Wolff/Phoenix	Wolff [®] base half-fat cream	609950	1
DR. WOLFF	Wolff [®] base half-fat cream	609950	1
Dr. August Wolff	Wolff [®] base half-fat cream	609950	1
Dr Wolff	Wolff [®] base half-fat cream	610991	1
Dr. Wolff/Kehr	Wolff [®] base half-fat cream	702500	3
Phönix WE: 22.09.17 EK...	Wolff [®] base half-fat cream	702500	1
Dr. Wolff/Kehr	Wolff [®] base half-fat cream	609950	2
Dr. Wolff/Phoenix	Wolff [®] base half-fat cream	702500	1
Sanacorp	Wolff [®] base half-fat cream	702500	1
Chiesi GmbH	Wolff [®] base half-fat cream	703031	1
Dr. Wolff/Sanacorp	Wolff [®] base half-fat cream	704050	3
Wolff	Wolff [®] base half-fat cream	704050	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Dr. August Wolff	Wolff [®] base half-fat cream	704050	1
Dr. August Wolff GmbH &...	Wolff [®] base half-fat cream		1
Dr Wolff	Wolff [®] base half-fat cream	704220	1
Dr. Wolff/Kehr	Wolff [®] base half-fat cream	707370	2
Wolff	Wolff [®] base half-fat cream	801140	1
Dr. Wolff GmbH/Phönix	Wolff [®] base half-fat cream	711250	1
Dr. Wolff	Wolff [®] base half-fat cream	10112014F	1
Dr Wolff/ Gehe	Wolff [®] base half-fat cream	234791	1
Caelo	Wolff [®] base half-fat cream	235200	1
Dr. Wolff	Wolff [®] base half-fat cream	331551	2
Wolff/Gehe	Wolff [®] base half-fat cream	236420	1
wolff/AHD	Wolff [®] base half-fat cream	332040	10
Wolff/AHD	Wolff [®] base half-fat cream	332040	5
Euro OTC	Wolff [®] base half-fat cream	332040	4
Caelo	Wolff [®] base half-fat cream	332040	6
Wolff	Wolff [®] base half-fat cream	332040	5
Phönix Gotha	Wolff [®] base half-fat cream	332040	2
Dr.Wolff/Sanacorp	Wolff [®] base half-fat cream	333620	1
Wolff GmbH/Sanacorp	Wolff [®] base half-fat cream	333620	3
Dr Wolff/ Gehe	Wolff [®] base half-fat cream	333930	1
Phoenix/Wolff	Wolff [®] base half-fat cream	333930	1
Dr. Wolff	Wolff [®] base half-fat cream	333930	3
Dr.Wolff/Noweda	Wolff [®] base half-fat cream	333930	2
Dr.Wolff	Wolff [®] base half-fat cream	339140	1
DR. WOLFF	Wolff [®] base half-fat cream	337260	1
VDL 05.05.14 Fa.Wolff	Wolff [®] base half-fat cream	339140	1
Wolff GmbH/Sanacorp	Wolff [®] base half-fat cream	339140	1
Dr. Wolff	Wolff [®] base half-fat cream	337260	3
Dr Wolff/ Sanacorp	Wolff [®] base half-fat cream	339140	1
Dr.Wolff	Wolff [®] base half-fat cream	337260	1
Cealo/Phönix	Wolff [®] base half-fat cream	339140	1
Anzag/ Dr.Wolff	Wolff [®] base half-fat cream	339140	1
VDL 28.04.2014;Preis:45...	Wolff [®] base half-fat cream		1
Wolff/Alliance	Wolff [®] base half-fat cream	430880	1
Wolf/Phönix	Wolff [®] base half-fat cream	430880	2
Dr. Wolff/Noweda	Wolff [®] base half-fat cream	430880	4
Hedinger/Noweda	Wolff [®] base half-fat cream	430880	2
VDL;2.7.14	Wolff [®] base half-fat cream	430880	3
Wolff	Wolff [®] base half-fat cream	430880	1
Wolff/ Phönix	Wolff [®] base half-fat cream	430880	2
Caelo	Wolff [®] base half-fat cream	431050	1
Dr.August Wolff	Wolff [®] base half-fat cream	430880	1
Dr. Wolff/ Gehe	Wolff [®] base half-fat cream	430880	1
Euro/Gehe	Wolff [®] base half-fat cream	431050	1
Gehe	Wolff [®] base half-fat cream	431050	1
Dr.Wolff	Wolff [®] base half-fat cream	431240	2
wolff	Wolff [®] base half-fat cream	431050	1
Caelo	Wolff [®] base half-fat cream	431240	2
Wolff	Wolff [®] base half-fat cream	431240	2
wolff/AHD	Wolff [®] base half-fat cream	431240	1
Dr. Wolff	Wolff [®] base half-fat cream	431240	1
wolff	Wolff [®] base half-fat cream	431240	1
Dr. Wolff /Anzag	Wolff [®] base half-fat cream	432611	1
Dr.wolff	Wolff [®] base half-fat cream	432611	1
Dr.Wolff GmbH & Co.KG A...	Wolff [®] base half-fat cream	432611	1
Dr.Wolff/NOW	Wolff [®] base half-fat cream	432611	1
Caelo	Wolff [®] base half-fat cream	434360	3
Dr. Wolff/ Sanacorp	Wolff [®] base half-fat cream	434360	1
Phoenix	Wolff [®] base half-fat cream	434360	1
Dr.Wolff/NOW	Wolff [®] base half-fat cream	434581	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Dr.Wolff GmbH & Co.KG A...	Wolff [®] base half-fat cream	434581	1
Dr.August Wolff	Wolff [®] base half-fat cream	38800	1
Phönix Gotha	Wolff [®] base half-fat cream	134990	1
Wolff GmbH	Wolff [®] base half-fat cream	230870	1
Löwen City-Apotheke Bar...	Wolff [®] base half-fat cream	230200	1
Spirig/Anzag	Wolff [®] base half-fat cream	234791	1
Dr. Wolff /Anzag	Wolff [®] base half-fat cream	232593	1
Dr.Wolff/Noweda	Wolff [®] base half-fat cream	234791	2
Sanacorp	Wolff [®] base half-fat cream	235200	2
Dr.Wolff/ Anzag	Wolff [®] base half-fat cream	235200	1
DR. August Wolff GmbH	Wolff [®] base half-fat cream	235200	1
Fagron	Wolff [®] base half-fat cream	236320	3
Dr.Wolff	Wolff [®] base half-fat cream	236320	1
Dr.Wolff	Wolff [®] base half-fat cream	236420	1
Dr. Wolff/Anzag	Wolff [®] base half-fat cream	238490	2
Fagron	Wolff [®] base half-fat cream	238560	2
Dr. Wolff	Wolff [®] base half-fat cream	238560	5
Caelo	Wolff [®] base half-fat cream	238560	1
Dr.Wolff	Wolff [®] base half-fat cream	238560	2
Wolff GmbH	Wolff [®] base half-fat cream	238560	1
Wolff/Sanacorp	Wolff [®] base half-fat cream	238560	1
Phönix	Wolff [®] base half-fat cream	238560	1
Dr.Wolff/Gehe	Wolff [®] base half-fat cream	238660	1
Dr.Wolff/Noweda	Wolff [®] base half-fat cream	238660	1
Wolff	Wolff [®] base half-fat cream	238660	1
Wolff/AHD	Wolff [®] base half-fat cream	331551	8
Dr. Wolf	Wolff [®] base half-fat cream	331551	1
Gehe	Wolff [®] base half-fat cream	331551	2
Dr.Wolff	Wolff [®] base half-fat cream	332040	1
Dr. Wolff/Sanacorp	Wolff [®] base half-fat cream	333620	1
Phönix	Wolff [®] base half-fat cream	333620	2
Dr.August Wolff	Wolff [®] base half-fat cream	333620	1
Dr.Wolff	Wolff [®] base half-fat cream	333620	5
Caelo	Wolff [®] base half-fat cream	333620	3
Dr. Wolff / Noweda	Wolff [®] base half-fat cream	333930	1
Dr.Wolff	Wolff [®] base half-fat cream	333930	2
Caelo	Wolff [®] base half-fat cream	333930	1
Dr. August Wolff/Sanaco...	Wolff [®] base half-fat cream	337260	1
Dr. Wolff	Wolff [®] base half-fat cream	432611	1

- 81 502 spectra from 501 *Apo-Ident* customers from a total of 24 079 batches from a further 323 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Hydrophilic ointments* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Hydrophilic ointments* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	2290	0	33 254
Type B	0	14 173	0	79 111
Type C	0	9602	248	81 350

The substance/substance group *Hydrophilic ointments* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9702 %)	100.0000 % (> 99.7380 %)
Type B	100.0000 % (> 99.9764 %)	100.0000 % (> 99.9577 %)
Type C	100.0000 % (> 99.6730 %)	97.4822 % (> 97.4518 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30631	30837	9.06	34.02
30758	31907	6.73	35.04
30769	30769	0.00	25.80
30837	30837	0.00	32.21
30880	31741	6.47	29.23
30895	31054	6.88	32.68
30931	31478	10.96	30.80
31002	31907	3.15	34.49
31008	31008	0.00	20.82
31039	31039	0.00	28.32
31040	31503	6.72	27.61
31042	31008	8.98	22.98
31054	31054	0.00	30.75
31086	30837	13.77	30.62
31188	31188	0.00	25.96
31231	31231	0.00	34.09
31241	31741	3.75	36.01
31440	31440	0.00	27.04
31450	32123	5.93	24.93
31470	31231	3.11	35.11
31477	31188	3.28	24.79
31478	31478	0.00	25.96
31487	30837	17.93	29.84
31500	31500	0.00	27.69
31503	31503	0.00	30.55

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31537	31537	0.00	34.23
31565	31996	8.86	29.27
31568	30769	5.04	26.04
31591	31996	6.42	25.07
31593	31996	3.34	26.22
31599	31188	2.29	25.84
31600	31039	7.03	30.00
31644	31799	3.44	36.13
31661	31537	8.77	34.26
31741	31741	0.00	34.34
31799	31799	0.00	35.37
31906	31906	0.00	32.35
31907	31907	0.00	34.11
31962	31996	13.08	34.01
31996	31996	0.00	26.58
32022	32022	0.00	28.71
32039	32039	0.00	21.42
32052	32052	0.00	28.50
32123	32123	0.00	27.06
32132	32132	0.00	30.25
32142	32142	0.00	20.24
32303	31799	2.53	33.61
32346	31008	5.29	21.59
32435	32435	0.00	27.11
32620	32620	0.00	36.53
32621	32621	0.00	21.06
32882	32882	0.00	38.68
33295	31503	6.36	26.67
33354	33354	0.00	36.18
33595	33595	0.00	35.17
33641	31906	12.85	30.36
33964	31799	8.57	37.42

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Hydroxyethyl salicylate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30638-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Hydroxyethyl salicylate; Hydroxyethyl salicylicum; Hydroxyethyl salicylas

Special notes

When selecting the *Hydroxyethyl salicylate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Hydroxyethyl salicylate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Hydroxyethyl sal...	12049001	30638	40	1402192
Caelo	Hydroxyethyl sal...	13285709	31339	40	20130918

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Hydroxyethyl salicylate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Hydroxyethyl salicylate*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Hydroxyethyl salicylate	16131002	32588	40
Caelo	Hydroxyethyl salicylate	16131009	33241	40
Caelo	Hydroxyethyl salicylate	172740	34226	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 3 spectra from 3 *Apo-Ident* customers from 3 batches from the substance/substance group *Hydroxyethyl salicylate*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Hydroxyethyl salicylate	13285702	1
Caelo	Hydroxyethyl salicylate	15414906	1
Caelo	Hydroxyethyl salicylate	2091302	1

- 91 349 spectra from 501 *Apo-Ident* customers from a total of 27 006 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Hydroxyethyl salicylate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Hydroxyethyl salicylate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	120	0	93 754
Type C	0	3	0	91 349

The substance/substance group *Hydroxyethyl salicylate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6615 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30638	30638	0.00	220.96
31339	30638	8.65	225.20

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Isopropyl alcohol 100%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30759-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Isopropyl alcohol 100%; Alcohol isopropylicus 100%

Special notes

When selecting the *Isopropyl alcohol 100%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Isopropyl alcohol 100%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Isopropyl alcoho...	12321002	30759	40	20121102
Caelo	Isopropyl alcoho...	13249425	31232	40	20130726
Caelo	Isopropyl alcoho...	14051031	31386	40	20140214
apomix	Isopropyl alcoho...	4103112	31517	40	1505476
Hedinger	Isopropyl alcoho...	070835B	32817	40	20160712

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Isopropyl alcohol 100%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 344 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Isopropyl alcohol 100%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
zu erfassen	Isopropyl alcohol 100%	zu erfassen	30101	40
apomix	Isopropyl alcohol 100%	4103112	31517 [†]	20
Hedinger	Isopropyl alcohol 100%	070922B	33256	40
Caelo	Isopropyl alcohol 100%	17300401	33938	40

- 93 734 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1734 spectra from 328 *Apo-Ident* customers from 775 batches from the substance/substance group *Isopropyl alcohol 100%*.
- These include spectra of independent samples from 770 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Hedinger / Noweda	Isopropyl alcohol 100%	701011B	1
Chem.Pharm.Labor.Dr.See...	Isopropyl alcohol 100%	70577B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70598B	2
Hedinger	Isopropyl alcohol 100%	70630B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70734B	2
Hedinger	Isopropyl alcohol 100%	70730B	2
Hedinger/phoenix	Isopropyl alcohol 100%	707488	1
Hedinger/Gehe	Isopropyl alcohol 100%	70737B	1
Hedinger	Isopropyl alcohol 100%	70788B	1
Hedinger / Nowed	Isopropyl alcohol 100%	70730B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70703B	1
phönix	Isopropyl alcohol 100%	70799B	1
Hedinger	Isopropyl alcohol 100%	70837B	1
Hedinger	Isopropyl alcohol 100%	70879B	2
GEHE	Isopropyl alcohol 100%	70882B	1
Phönix	Isopropyl alcohol 100%	70862B	1
Hedinger/VDL	Isopropyl alcohol 100%	70884B	1
Allianz Healthcare	Isopropyl alcohol 100%	70896b	1
GEHE	Isopropyl alcohol 100%	70896B	1
AHD	Isopropyl alcohol 100%	70898B	1
Sanacorp	Isopropyl alcohol 100%	70918B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70922B	2
Hedinger	Isopropyl alcohol 100%	70918B	1
Hedinger/noweda	Isopropyl alcohol 100%	70957B	1
Noweda	Isopropyl alcohol 100%	70993B	1
Fischar	Isopropyl alcohol 100%	1005101	1
Gehe/Fischar	Isopropyl alcohol 100%	1005101	2
Fischar	Isopropyl alcohol 100%	1007011	1
Caelo	Isopropyl alcohol 100%	13336927	1
Caelo	Isopropyl alcohol 100%	13397327	2
Caelo	Isopropyl alcohol 100%	14219212	3
Caelo	Isopropyl alcohol 100%	14051022	1
Caelo	Isopropyl alcohol 100%	14278904	3
Caelo	Isopropyl alcohol 100%	14219216	1
Caelo	Isopropyl alcohol 100%	14278914	2
Caelo	Isopropyl alcohol 100%	14369312	2
Caelo	Isopropyl alcohol 100%	14369313	1
Caelo	Isopropyl alcohol 100%	1469323	1
Caelo	Isopropyl alcohol 100%	15300902	1
Caelo	Isopropyl alcohol 100%	15398827	2
Ebert & Jacobi 28.7.16	Isopropyl alcohol 100%	16137317	1
Caelo	Isopropyl alcohol 100%	16137325	2
Caelo	Isopropyl alcohol 100%	16137329	1
Caelo	Isopropyl alcohol 100%	16234032	1
Caelo	Isopropyl alcohol 100%	16310622	3
Gehe	Isopropyl alcohol 100%	17044002	1
Caelo	Isopropyl alcohol 100%	17044030	2
Noweda	Isopropyl alcohol 100%	17044033	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	17044035	1
Caelo	Isopropyl alcohol 100%	17189302	1
Caelo	Isopropyl alcohol 100%	17189303	1
Caelo	Isopropyl alcohol 100%	17300412	1
Caelo	Isopropyl alcohol 100%	17300422	1
GEHE	Isopropyl alcohol 100%	17300433	1
Caelo	Isopropyl alcohol 100%	18065402	1
GEHE	Isopropyl alcohol 100%	18065406	1
fischar	Isopropyl alcohol 100%	3346271	1
Hedinger	Isopropyl alcohol 100%	4103014	1
Hedinger , Noweda	Isopropyl alcohol 100%	70490B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70501B	1
Otto Fischar	Isopropyl alcohol 100%	ohne	1
Noweda/Hedinger	Isopropyl alcohol 100%	70507B	9
Hedinger	Isopropyl alcohol 100%	705258	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70540B	1
Caelo	Isopropyl alcohol 100%	70541b	2
Hedinger/Noweda	Isopropyl alcohol 100%	70620B	1
Caelo	Isopropyl alcohol 100%	12179316	2
Caelo	Isopropyl alcohol 100%	12254703	1
Caelo	Isopropyl alcohol 100%	12031116	1
Caelo	Isopropyl alcohol 100%	12254727	2
Phönix	Isopropyl alcohol 100%	12321002	2
Sanacorp	Isopropyl alcohol 100%	12321028	1
Caelo	Isopropyl alcohol 100%	12321030	2
Caelo	Isopropyl alcohol 100%	13148004	1
Caelo	Isopropyl alcohol 100%	12379504	1
Caelo	Isopropyl alcohol 100%	13336913	1
Caelo	Isopropyl alcohol 100%	13397305	2
Caelo	Isopropyl alcohol 100%	30011401	1
Phönix	Isopropyl alcohol 100%	15300908	1
Caelo	Isopropyl alcohol 100%		5
Fischar	Isopropyl alcohol 100%	4100111	1
Kassner & Sasse	Isopropyl alcohol 100%	383316	2
Hedinger	Isopropyl alcohol 100%	70562B	3
Hedinger	Isopropyl alcohol 100%	70564B	2
Hedinger/Phönix	Isopropyl alcohol 100%	70565B	3
Hedinger	Isopropyl alcohol 100%	70593B	3
Hedinger	Isopropyl alcohol 100%	70598B	9
Hedinger/ N/ PZN:07474942	Isopropyl alcohol 100%	70598B	1
Gutenberg-Apotheke	Isopropyl alcohol 100%	70607B	1
Hedinger	Isopropyl alcohol 100%	70636B	2
Hedinger/Phönix	Isopropyl alcohol 100%	70643B	1
Aug. hedinger/N/PZN:074...	Isopropyl alcohol 100%	70659B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70671b	1
Hedinger Gehe	Isopropyl alcohol 100%	70673B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70675B	1
Hedinger/Phönix	Isopropyl alcohol 100%	706808	1
Hedinger/Anzag	Isopropyl alcohol 100%	70685B	1
Hedinger Jenne	Isopropyl alcohol 100%	706938	1
Hedinger/Gehe	Isopropyl alcohol 100%	70703B	1
Hedinger/Phönix	Isopropyl alcohol 100%	70715B	2
Hedinger	Isopropyl alcohol 100%	70715b	1
Hedinger vdL	Isopropyl alcohol 100%	70715B	1
Hedinger Jenne	Isopropyl alcohol 100%	70720B	1
Hedinger/Phoenix	Isopropyl alcohol 100%	70799B	1
Hedinger vdL	Isopropyl alcohol 100%	70748B	1
Gehe/Hedinger	Isopropyl alcohol 100%	70802B	1
Hedinger/Phoenix	Isopropyl alcohol 100%	70734B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70763B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Phoenix	Isopropyl alcohol 100%	70794B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70818B	1
Hedinger vdL	Isopropyl alcohol 100%	70853B	1
sanacorp, 29,30EUR, 15.09...	Isopropyl alcohol 100%	70835B	1
HedingerAlliance	Isopropyl alcohol 100%	70866B	2
HedingerAlliance	Isopropyl alcohol 100%	70881B	2
HedingerAlliance	Isopropyl alcohol 100%	70882B	1
Hedinger vdL	Isopropyl alcohol 100%	70889B	1
HEDINGER / NOWEDA	Isopropyl alcohol 100%	70898B	1
VDL	Isopropyl alcohol 100%	14219213	1
Sanacorp, 18.03.17, 29,...	Isopropyl alcohol 100%	70906B	1
Caelo	Isopropyl alcohol 100%	14051016	1
Caelo	Isopropyl alcohol 100%	14219227	3
Caelo	Isopropyl alcohol 100%	14051012	3
Caelo	Isopropyl alcohol 100%	14278905	2
Caelo	Isopropyl alcohol 100%	14369302	5
Caelo	Isopropyl alcohol 100%	14278923	4
Caelo	Isopropyl alcohol 100%	14278926	4
Caelo	Isopropyl alcohol 100%	14219233	1
Caelo	Isopropyl alcohol 100%	15057501	2
Caelo	Isopropyl alcohol 100%	15057502	2
Caelo	Isopropyl alcohol 100%	14369311	1
Caelo	Isopropyl alcohol 100%	14369325	1
Caelo	Isopropyl alcohol 100%	15171304	3
Caelo	Isopropyl alcohol 100%	15171329	2
Caelo	Isopropyl alcohol 100%	15398828	1
Caelo	Isopropyl alcohol 100%	16137302	2
Caelo	Isopropyl alcohol 100%	16039911	2
Caelo	Isopropyl alcohol 100%	16039923	2
Caelo	Isopropyl alcohol 100%	16137307	1
Caelo	Isopropyl alcohol 100%	16137317	1
Caelo	Isopropyl alcohol 100%	162334019	1
Caelo	Isopropyl alcohol 100%	16234009	1
Caelo	Isopropyl alcohol 100%	16234013	1
Caelo	Isopropyl alcohol 100%	17044027	1
Otto Fischar GmbH	Isopropyl alcohol 100%	4104072	1
Noweda	Isopropyl alcohol 100%	70409B	1
Gehe	Isopropyl alcohol 100%	70448B	1
GEHE	Isopropyl alcohol 100%	70461B	1
Anzag	Isopropyl alcohol 100%	70471B	1
Anzag	Isopropyl alcohol 100%	70459B	1
Caelo	Isopropyl alcohol 100%	70481B	2
Hedinger/ Anzag	Isopropyl alcohol 100%	70492B	1
Herbasinica	Isopropyl alcohol 100%	70478b	1
Jenne, Hedinger	Isopropyl alcohol 100%	70501B	1
Anzag	Isopropyl alcohol 100%	70495B	1
Jenne, Hedinger	Isopropyl alcohol 100%	70519B	1
Hedinger/Gehe	Isopropyl alcohol 100%	704988	1
Hedinger/Phönix	Isopropyl alcohol 100%	70520B	1
Anzag/Hedinger	Isopropyl alcohol 100%	70495B	1
HEDINGER/N/1.06.13/PZN...	Isopropyl alcohol 100%	70527b	1
SHELL	Isopropyl alcohol 100%	70540b	1
Hedinger/AHD	Isopropyl alcohol 100%	705408	1
Hedinger Jenne	Isopropyl alcohol 100%	70545B	1
Hedinger/N/PZN:07474942	Isopropyl alcohol 100%	70559B	1
Hedinger/AHD	Isopropyl alcohol 100%	70527B	1
Hedinger Gehe	Isopropyl alcohol 100%	70562B	1
Hedinger/AHD	Isopropyl alcohol 100%	70574B	1
Hedinger/Gehe	Isopropyl alcohol 100%	7087B	1
Caelo	Isopropyl alcohol 100%	12379526	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	13050913	2
Sanacorp	Isopropyl alcohol 100%	13050904	1
Caelo	Isopropyl alcohol 100%	132494	1
Caelo	Isopropyl alcohol 100%	13148001	2
Caelo	Isopropyl alcohol 100%	13249406	3
Sanacorp	Isopropyl alcohol 100%	13249410	1
Fischar	Isopropyl alcohol 100%	4102084	2
Hedinger	Isopropyl alcohol 100%	2072015A	1
Hedinger	Isopropyl alcohol 100%	2022015D	1
Caelo	Isopropyl alcohol 100%	4061401	1
Caelo	Isopropyl alcohol 100%	4031405	1
Caelo	Isopropyl alcohol 100%	5061415	1
Caelo	Isopropyl alcohol 100%	53380	1
Caelo	Isopropyl alcohol 100%	6061402	1
Caelo	Isopropyl alcohol 100%	6061403	1
Caelo	Isopropyl alcohol 100%	6061404	1
Hedinger/Noweda	Isopropyl alcohol 100%	701006B	1
Hedinger	Isopropyl alcohol 100%	7010318	1
hedinger	Isopropyl alcohol 100%	701011b	1
Hedinger/L10	Isopropyl alcohol 100%	7010158	1
Hedinger / Phoenix	Isopropyl alcohol 100%	70536B	1
Caelo	Isopropyl alcohol 100%	7051401	1
Hedinger/Gehe	Isopropyl alcohol 100%	70562B	1
Hedinger , Stuttgart	Isopropyl alcohol 100%	70545B	1
Hedinger	Isopropyl alcohol 100%	70572B	5
Caelo	Isopropyl alcohol 100%	705538	2
Hedinger/Phönix	Isopropyl alcohol 100%	70536B	1
Gehe am 07.03.14	Isopropyl alcohol 100%	70573B	1
Hedinger(Gehe)	Isopropyl alcohol 100%	70573B	1
Sanacorp / Hedinger	Isopropyl alcohol 100%	70573B	1
Caelo	Isopropyl alcohol 100%	70573B	1
Hedinger/Nowenda	Isopropyl alcohol 100%	70577B	2
Shell Chemicals	Isopropyl alcohol 100%	70577B	1
Gehe am 17.12.2013	Isopropyl alcohol 100%	70578B	1
Alliance Healthcare/Hed...	Isopropyl alcohol 100%	70578B	1
Hedinger / Gehe	Isopropyl alcohol 100%	70581B	1
Hedinger/Phoenix	Isopropyl alcohol 100%	70578B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70581B	1
Hedinger	Isopropyl alcohol 100%	70584B	3
Caelo	Isopropyl alcohol 100%	70584B	1
Ichthyl/Noweda	Isopropyl alcohol 100%	70587B	1
Hedinger / Apotheke E.B...	Isopropyl alcohol 100%	70587B	1
Hedinger	Isopropyl alcohol 100%	70592B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70587B	1
Hedinger	Isopropyl alcohol 100%	70592B	1
Anzag/Hedinger	Isopropyl alcohol 100%	70592B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70593B	1
Sanacorp	Isopropyl alcohol 100%	705988	1
Hedinger/ Jenne	Isopropyl alcohol 100%	70595B	1
Caelo	Isopropyl alcohol 100%	705988	1
Anzag/Hedinger	Isopropyl alcohol 100%	70598B	4
Noweda / Hedinger	Isopropyl alcohol 100%	70598B	1
Caelo	Isopropyl alcohol 100%	70598B	4
Caelo	Isopropyl alcohol 100%	70598b	1
Isopropylalkohol GMP/Gehe	Isopropyl alcohol 100%	70598B	1
GEHE 08.05.2014	Isopropyl alcohol 100%	70598B	1
Hedinger/Phönix	Isopropyl alcohol 100%	70598B	1
Hedinger	Isopropyl alcohol 100%	706034	1
Hedinger	Isopropyl alcohol 100%	70603B	6
Gehe/Hedinger	Isopropyl alcohol 100%	70603B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Noweda	Isopropyl alcohol 100%	70603B	3
Hedinger/Sanacorp6.5.2014	Isopropyl alcohol 100%	70603B	1
	Isopropyl alcohol 100%	70603B	1
Hedinger/Alliance	Isopropyl alcohol 100%	70603B	3
Hedinger/Gehe	Isopropyl alcohol 100%	70603B	1
Hedinger	Isopropyl alcohol 100%	70607B	3
Hedinger/Holdermann	Isopropyl alcohol 100%	70607B	1
Hedinger1	Isopropyl alcohol 100%	70607B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70607B	1
Rondell Apotheke	Isopropyl alcohol 100%	70607B	2
Hedinger/Phönix	Isopropyl alcohol 100%	70607B	2
Gehe - Hedinger	Isopropyl alcohol 100%	70607B	1
Phoenix 22.10.14	Isopropyl alcohol 100%	70607B	1
Hedinger/ Jenne	Isopropyl alcohol 100%	70608B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70608B	1
Hedinger	Isopropyl alcohol 100%	7060B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70614B	3
Alliance/ Hedinger	Isopropyl alcohol 100%	70610B	1
AHD	Isopropyl alcohol 100%	70614B	1
A. Hedinger / Kehr	Isopropyl alcohol 100%	70614B	1
HEDINGER / 20	Isopropyl alcohol 100%	70614B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70614B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70620B	1
Hedinger	Isopropyl alcohol 100%	70620B	6
Hedinger , Stuttgart	Isopropyl alcohol 100%	70620B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70620B	1
Hedinger/ Jenne	Isopropyl alcohol 100%	70620B	1
Hedinger	Isopropyl alcohol 100%	70623B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70627B	2
Hedinger/Sanacorp	Isopropyl alcohol 100%	70627B	3
Hedinger/ Noweda	Isopropyl alcohol 100%	70627B	2
Hedinger	Isopropyl alcohol 100%	70627B	6
Noweda/Hedinger	Isopropyl alcohol 100%	70627B	1
Hedinger/Phönix	Isopropyl alcohol 100%	70627B	1
Alliance Healthcare/Hed. . .	Isopropyl alcohol 100%	70627B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70627B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70627B	1
Noweda	Isopropyl alcohol 100%	70627B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70627B	3
Stumpf/Hedinger	Isopropyl alcohol 100%	70627B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70638B	1
Euro OTC	Isopropyl alcohol 100%	70638	1
Hedinger/Noweda	Isopropyl alcohol 100%	706388	1
Hedinger	Isopropyl alcohol 100%	70638B	3
Hedinger/Geilenkirchen	Isopropyl alcohol 100%	706388	1
Hedinger/ GEHE	Isopropyl alcohol 100%	706388	1
Noweda / Hedinger	Isopropyl alcohol 100%	70638B	1
Hedinger/ Anzag	Isopropyl alcohol 100%	706388	1
Hedinger/Jenne 31.10.2014	Isopropyl alcohol 100%	706388	1
Noweda	Isopropyl alcohol 100%	70638B	1
Hedinger	Isopropyl alcohol 100%	70643B	9
Hedinger/Gehe	Isopropyl alcohol 100%	70638B	1
HEDINGER / 20	Isopropyl alcohol 100%	70638B	1
Jenne/Hedinger	Isopropyl alcohol 100%	70643B	1
Caelo	Isopropyl alcohol 100%	70643B	2
Hedinger / Apotheke E.B. . .	Isopropyl alcohol 100%	70643B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70643B	1
AHD	Isopropyl alcohol 100%	70643B	1
hedinger	Isopropyl alcohol 100%	70643B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70650B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Sanacorp	Isopropyl alcohol 100%	70650B	1
Hedinger	Isopropyl alcohol 100%	70650B	6
Caelo	Isopropyl alcohol 100%	70650B	1
Hedinger	Isopropyl alcohol 100%	70656B	1
Caelo	Isopropyl alcohol 100%	70656B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70656B	1
Hedinger	Isopropyl alcohol 100%	70659B	4
Noweda/Hedinger	Isopropyl alcohol 100%	70659B	1
A. Hedinger / Noweda	Isopropyl alcohol 100%	70659B	2
Hedinger/Sanacorp	Isopropyl alcohol 100%	70659B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70659B	2
Hedinger / Noweda	Isopropyl alcohol 100%	70659B	1
A.Hedinger / Kehr	Isopropyl alcohol 100%	70659B	1
Hedinger/Phoenix	Isopropyl alcohol 100%	70659B	2
Hedinger	Isopropyl alcohol 100%	70663B	3
Hedinger / Noweda	Isopropyl alcohol 100%	70663B	1
Noweda/Henry Lamotte	Isopropyl alcohol 100%	70663B	1
hedinger	Isopropyl alcohol 100%	70663B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70663B	1
Noweda/ Hedinger GmbH	Isopropyl alcohol 100%	70663B	1
AHD /Hedinger	Isopropyl alcohol 100%	70666B	1
Hedinger	Isopropyl alcohol 100%	70666B	1
Caelo	Isopropyl alcohol 100%	70663B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70666B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70668B	1
Hedinger	Isopropyl alcohol 100%	70668B	1
Hedinger	Isopropyl alcohol 100%	70671B	2
Caelo	Isopropyl alcohol 100%	70668B	1
Hedinger	Isopropyl alcohol 100%	70675B	3
Hedinger/Kehr	Isopropyl alcohol 100%	70675B	1
Stumpf/Hedinger	Isopropyl alcohol 100%	70673B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70675B	1
Gehe-Hedinger	Isopropyl alcohol 100%	70675B	1
Hediger/ Noweda	Isopropyl alcohol 100%	70675B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70675B	1
hedinger	Isopropyl alcohol 100%	70675B	1
Hedinger/N	Isopropyl alcohol 100%	70675B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70680B	1
Alliance Healthcare/Hed...	Isopropyl alcohol 100%	70680B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70675B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70685	1
Hedinger	Isopropyl alcohol 100%	70680B	1
Hedinger/Phönix	Isopropyl alcohol 100%	70680B	1
Caelo	Isopropyl alcohol 100%	70685B	2
Hedinger	Isopropyl alcohol 100%	70685B	6
Hedinger/Sanacorp	Isopropyl alcohol 100%	70685B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70685B	3
Phoenix goettin	Isopropyl alcohol 100%	70685B	1
Fagron	Isopropyl alcohol 100%	70685B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70685B	1
Gehe	Isopropyl alcohol 100%	70685B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70690B	2
Hedinger	Isopropyl alcohol 100%	70690B	3
Caelo	Isopropyl alcohol 100%	70690B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70693B	2
Hedinger	Isopropyl alcohol 100%	70693B	2
Hedinger/ Noweda	Isopropyl alcohol 100%	70693B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70693B	1
Caelo	Isopropyl alcohol 100%	70693B	1
Hedinger / AHD 21.4.15	Isopropyl alcohol 100%	70693B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
August Hedinger/Phoenix	Isopropyl alcohol 100%	70695B	1
Hedinger	Isopropyl alcohol 100%	70695B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70696B	1
Hedinger	Isopropyl alcohol 100%	70701B	2
Hedinger/Sanacorp	Isopropyl alcohol 100%	70696B	1
Sanacorp	Isopropyl alcohol 100%	70701B	1
Eigen	Isopropyl alcohol 100%	70701B	1
Hedinger	Isopropyl alcohol 100%	70703B	3
Hedinger/ Kehr	Isopropyl alcohol 100%	70703B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70703B	1
Aug.Hedinger / Phoenix	Isopropyl alcohol 100%	70703B	1
Shell Chemicals	Isopropyl alcohol 100%	70703B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70703B	1
Caelo	Isopropyl alcohol 100%	70707B	1
Hedinger	Isopropyl alcohol 100%	70707B	6
Hedinger	Isopropyl alcohol 100%	70704B	1
Hedinger	Isopropyl alcohol 100%	70708B	2
Hedinger/ Gehe	Isopropyl alcohol 100%	70708B	2
A. Hedinger / Kehr	Isopropyl alcohol 100%	707158	1
Hedinger	Isopropyl alcohol 100%	707158	1
Hedinger	Isopropyl alcohol 100%	70715B	5
Hedinger/Gehe	Isopropyl alcohol 100%	70715B	1
Hedinger	Isopropyl alcohol 100%	70720B	2
Audor Pharma	Isopropyl alcohol 100%	70720B	2
Hedinger/Sanacorp	Isopropyl alcohol 100%	70720B	1
Anzag	Isopropyl alcohol 100%	70720B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70720B	1
Hedinger	Isopropyl alcohol 100%	70722B	3
Noweda/Hedinger	Isopropyl alcohol 100%	70722B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70723B	1
Hedinger	Isopropyl alcohol 100%	70723B	1
Hedinger/Phönix	Isopropyl alcohol 100%	70723B	1
HEDINGER / 20	Isopropyl alcohol 100%	70726B	1
Noweda/Hedin	Isopropyl alcohol 100%	70726B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70726B	1
Caelo	Isopropyl alcohol 100%	70726B	1
Hedinger/Gehe	Isopropyl alcohol 100%	707308	1
Noweda/Hedinger	Isopropyl alcohol 100%	70730B	1
Hedinger/Phoenix	Isopropyl alcohol 100%	70730B	4
Allianc. Healthc. / Hed...	Isopropyl alcohol 100%	70730B	1
Heinger/G	Isopropyl alcohol 100%	70730B	1
Hedinger GmbH	Isopropyl alcohol 100%	70729B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70730B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70734B	1
Phönix 13.10.15 ...	Isopropyl alcohol 100%	70734B	1
Hedinger / Noweda Münster	Isopropyl alcohol 100%	70734B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70734B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70736B	1
Hediger/Sanacorp	Isopropyl alcohol 100%	70736B	1
Hedinger	Isopropyl alcohol 100%	70736B	2
Hedinger/Anzag	Isopropyl alcohol 100%	70736B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70737B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70737B	2
Hedinger/AHD	Isopropyl alcohol 100%	70737B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70737B	1
Hedinger / Sanacorp	Isopropyl alcohol 100%	70737B	1
Hedinger	Isopropyl alcohol 100%	70737B	2
Hedinger Sanacorp	Isopropyl alcohol 100%	70737B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70737B	1
Hedinger/sanacorp	Isopropyl alcohol 100%	70739B	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	70739B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70739B	1
Hedinger/Kehr	Isopropyl alcohol 100%	7073B	1
Hedinger	Isopropyl alcohol 100%	70739B	1
Noweda / Hedinger	Isopropyl alcohol 100%	70742B	2
Hedinger	Isopropyl alcohol 100%	70742B	2
Hedinger / Gehe	Isopropyl alcohol 100%	70742B	1
Hedinger	Isopropyl alcohol 100%	70748B	4
Heidinger	Isopropyl alcohol 100%	70748B	1
Caelo	Isopropyl alcohol 100%	70748B	2
Sanacorp/Hedinger WE:15...	Isopropyl alcohol 100%	70748B	1
Hedinger / Sana	Isopropyl alcohol 100%	70748B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70748B	2
Hedinger / Noweda	Isopropyl alcohol 100%	70748B	1
Hedinger/Phönix	Isopropyl alcohol 100%	70748B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70748b	1
Hedinger	Isopropyl alcohol 100%	70755B	2
Hedinger/Gehe	Isopropyl alcohol 100%	70755B	2
Noweda	Isopropyl alcohol 100%	70755B	1
Caelo	Isopropyl alcohol 100%	70755B	2
Hedinger / Kehr	Isopropyl alcohol 100%	70755B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70763B	2
Hedinger / Noweda	Isopropyl alcohol 100%	70762B	1
AHD/Hedinger	Isopropyl alcohol 100%	70763B	1
Hedinger	Isopropyl alcohol 100%	70763B	1
Allianc. Health. / Hedi...	Isopropyl alcohol 100%	70762B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70767B	5
Hedinger	Isopropyl alcohol 100%	70767B	3
Bombastus	Isopropyl alcohol 100%	70767B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70767B	2
Hedinger/Gehe	Isopropyl alcohol 100%	70767B	1
Hedinger/Alliance Healt...	Isopropyl alcohol 100%	70768B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70768B	1
Hedinger	Isopropyl alcohol 100%	70768B	1
Hedinger / 30	Isopropyl alcohol 100%	70768B	1
Caelo	Isopropyl alcohol 100%	70768b	1
Hedinger/Phönix	Isopropyl alcohol 100%	70768B	1
Aug. Hedinger GmbH & Co...	Isopropyl alcohol 100%	70768B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70771B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70771B	1
Hedinger	Isopropyl alcohol 100%	70771B	4
Hedinger/Sanacorp	Isopropyl alcohol 100%	70777B	1
Hedinger	Isopropyl alcohol 100%	70783B	3
Hedinger/Phönix	Isopropyl alcohol 100%	70783B	1
sittig	Isopropyl alcohol 100%	70788B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70783B	1
Hedinger/AHD	Isopropyl alcohol 100%	70788B	1
phönix	Isopropyl alcohol 100%	70783b	1
Hedinger/AHD	Isopropyl alcohol 100%	707888	1
Bartel/Gehe	Isopropyl alcohol 100%	70788B	2
Hedinger/Anzag	Isopropyl alcohol 100%	70788B	2
Hedinger/ AHD	Isopropyl alcohol 100%	70788B	1
Caelo	Isopropyl alcohol 100%	70788B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70788B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70788B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70788B	1
Hedinger / Sana	Isopropyl alcohol 100%	70799B	1
Hedinger	Isopropyl alcohol 100%	70794B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70799B	1
Hedinger	Isopropyl alcohol 100%	70799B	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Isopropyl alcohol 100%	70796B	1
Hedinger / Alliance Hea...	Isopropyl alcohol 100%	70796B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70796B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70799B	1
Hedinger	Isopropyl alcohol 100%	70802B	1
Hedinger / Gehe	Isopropyl alcohol 100%	70802B	1
Hedinger	Isopropyl alcohol 100%	70809B	6
Hedinger/Gehe	Isopropyl alcohol 100%	70804B	1
Sanacorp	Isopropyl alcohol 100%	70809B	3
Hedinger/Phönix	Isopropyl alcohol 100%	70809B	1
Phönix / Hedinger	Isopropyl alcohol 100%	70809B	1
Hedinger / Gehe	Isopropyl alcohol 100%	70809B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70809B	1
hedinger	Isopropyl alcohol 100%	70818B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70818B	4
Hedinger/Noweda	Isopropyl alcohol 100%	70818B	1
Hedinger	Isopropyl alcohol 100%	70818B	2
Hedinger/Gehe	Isopropyl alcohol 100%	70818B	3
Hedinger/N	Isopropyl alcohol 100%	70818B	1
Rondell Apotheke	Isopropyl alcohol 100%	70818B	1
Hedinger	Isopropyl alcohol 100%	70825B	4
Hedinger/Phönix	Isopropyl alcohol 100%	70825B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70825B	1
A. Hedinger / Kehr	Isopropyl alcohol 100%	70835B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70835B	1
Hedinger	Isopropyl alcohol 100%	70835B	4
Hedinger Sanacorp	Isopropyl alcohol 100%	70835B	1
Anzag/Hedinger	Isopropyl alcohol 100%	70835B	1
Gehe	Isopropyl alcohol 100%	70835B	1
Ichthyol/Phönix	Isopropyl alcohol 100%	70835B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70835B	2
Caelo	Isopropyl alcohol 100%	70835B	1
phönix	Isopropyl alcohol 100%	70837B	1
Allianc. Healthc. / Hed...	Isopropyl alcohol 100%	70837B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70853B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70853B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70853B	1
Hedinger	Isopropyl alcohol 100%	70853B	1
Hedinger	Isopropyl alcohol 100%	70858b	1
Hedinger	Isopropyl alcohol 100%	70858B	5
Hedinger	Isopropyl alcohol 100%	70855B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70858B	1
Hedinger Sanacorp	Isopropyl alcohol 100%	70862B	1
Hedinger	Isopropyl alcohol 100%	70862B	1
phönix	Isopropyl alcohol 100%	70862B	1
L 10/Hedinger	Isopropyl alcohol 100%	70866B	1
WE:Sanacorp 08.11.2016 ...	Isopropyl alcohol 100%	70870B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70875B	1
Phönix	Isopropyl alcohol 100%	70870B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70875B	2
Hedinger / Noweda	Isopropyl alcohol 100%	70870B	1
Caelo	Isopropyl alcohol 100%	70870B	1
Hedinger	Isopropyl alcohol 100%	70870B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70870B	1
Hedinger / Sana	Isopropyl alcohol 100%	70881B	1
Hedinger/Holdermann	Isopropyl alcohol 100%	70881B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70889B	1
Hedinger	Isopropyl alcohol 100%	70881B	1
Hedinger	Isopropyl alcohol 100%	70882B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70882b	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger / Noweda	Isopropyl alcohol 100%	70889B	1
Hedinger GmbH	Isopropyl alcohol 100%	70889B	1
Caelo	Isopropyl alcohol 100%	70889B	1
Hedinger	Isopropyl alcohol 100%	70889B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70892B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70892B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70889B	1
Hedinger	Isopropyl alcohol 100%	70898B	4
Hedinger/Holdermann	Isopropyl alcohol 100%	70898B	1
phönix	Isopropyl alcohol 100%	70898b	1
Hedinger	Isopropyl alcohol 100%	7089B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70908B	1
Hedinger	Isopropyl alcohol 100%	70908B	1
A. Hedinger/Kehr	Isopropyl alcohol 100%	709088	1
Hedinger/Noweda	Isopropyl alcohol 100%	70906B	1
Hedinger	Isopropyl alcohol 100%	70906B	1
Hedinger	Isopropyl alcohol 100%	70915B	1
Sanacorp EK:4,56 EUR -11	Isopropyl alcohol 100%	70923B	1
Hedinger	Isopropyl alcohol 100%	70922B	1
Hedinger	Isopropyl alcohol 100%	70931B	3
Hedinger Sanacorp	Isopropyl alcohol 100%	70931B(I)	1
Hedinger/Gehe	Isopropyl alcohol 100%	70931B	2
Hedinger Sanacorp	Isopropyl alcohol 100%	70931B(II)	2
Hedinger GmbH	Isopropyl alcohol 100%	70936B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70941B	4
Fagron	Isopropyl alcohol 100%	70936B	1
Hedinger	Isopropyl alcohol 100%	70936B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70943B	1
HEDINGER/Gehe	Isopropyl alcohol 100%	709478	1
Hedinger	Isopropyl alcohol 100%	70957B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70957B	1
Hedinger/ AHD	Isopropyl alcohol 100%	70947B	1
Hedinger	Isopropyl alcohol 100%	70973	1
Caelo	Isopropyl alcohol 100%	70959B	2
Hedinger	Isopropyl alcohol 100%	70959B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70973B	1
Hedinger	Isopropyl alcohol 100%	70977B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70977B	2
Phönix	Isopropyl alcohol 100%	70983B	1
Hedinger	Isopropyl alcohol 100%	70983B	2
Hedinger/ AHD	Isopropyl alcohol 100%	70991B	1
Hedinger	Isopropyl alcohol 100%	70991B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70994B	1
Phönix W.E:11.04.18 E...	Isopropyl alcohol 100%	70991B	1
Hedinger	Isopropyl alcohol 100%	70994B	2
Hedinger / Kehr	Isopropyl alcohol 100%	7098B	1
Hedinger / Phönix	Isopropyl alcohol 100%	70994B	1
Caelo	Isopropyl alcohol 100%	70994B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70994B	1
Bombastus	Isopropyl alcohol 100%	7474936	1
Hedinger, Noweda	Isopropyl alcohol 100%	7474907	1
Caelo	Isopropyl alcohol 100%	8082014C	1
Hedinger	Isopropyl alcohol 100%	7638B	1
Fischar	Isopropyl alcohol 100%		1
Caelo	Isopropyl alcohol 100%	10011502	1
Caelo	Isopropyl alcohol 100%	10061507	1
Caelo	Isopropyl alcohol 100%	10101403	1
Spinnrad	Isopropyl alcohol 100%	1114410092	1
Caelo	Isopropyl alcohol 100%	11091401	1
Caelo	Isopropyl alcohol 100%	12081505	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	12041402	1
Caelo	Isopropyl alcohol 100%	13121401	1
Caelo	Isopropyl alcohol 100%	13148006	1
Caelo	Isopropyl alcohol 100%	13336916	1
Caelo	Isopropyl alcohol 100%	13397303	1
Caelo	Isopropyl alcohol 100%	13397314	1
Caelo	Isopropyl alcohol 100%	13336921	1
Caelo	Isopropyl alcohol 100%	13397316	2
Caelo	Isopropyl alcohol 100%	13397312	2
Caelo	Isopropyl alcohol 100%	13397317	1
Caelo	Isopropyl alcohol 100%	13397318	2
Caelo	Isopropyl alcohol 100%	13397330	6
Euro OTC	Isopropyl alcohol 100%	13397324	1
Großhandel	Isopropyl alcohol 100%	13397331	1
Caelo	Isopropyl alcohol 100%	13397332	2
Caelo	Isopropyl alcohol 100%	13407304	1
Caelo	Isopropyl alcohol 100%	14031501	1
Caelo	Isopropyl alcohol 100%	14051008	1
Caelo	Isopropyl alcohol 100%	14051009	2
Sanacorp	Isopropyl alcohol 100%	14051017	1
Caelo	Isopropyl alcohol 100%	14051013	1
Caelo	Isopropyl alcohol 100%	14051017	1
Caelo	Isopropyl alcohol 100%	14051015	1
VDL	Isopropyl alcohol 100%	14051020	1
Caelo	Isopropyl alcohol 100%	14051018	1
Caelo	Isopropyl alcohol 100%	14051021	1
Caelo	Isopropyl alcohol 100%	14051031	3
Caelo	Isopropyl alcohol 100%	14051026	1
Noweda	Isopropyl alcohol 100%	14051027	1
Caelo	Isopropyl alcohol 100%	14051028	1
Caelo	Isopropyl alcohol 100%	140561033	1
Caelo	Isopropyl alcohol 100%	14051034	1
Sanacorp	Isopropyl alcohol 100%	14051035	1
Caelo	Isopropyl alcohol 100%	14051036	1
Caelo	Isopropyl alcohol 100%	141/AKD16	1
Caelo	Isopropyl alcohol 100%	14134703	1
Caelo	Isopropyl alcohol 100%	14134706	1
Caelo	Isopropyl alcohol 100%	14134716	3
Bergchemie	Isopropyl alcohol 100%	14134711	1
Caelo	Isopropyl alcohol 100%	14134705	2
Sanacorp	Isopropyl alcohol 100%	14134709	1
Caelo	Isopropyl alcohol 100%	14134710	1
Caelo	Isopropyl alcohol 100%	14134719	1
Fagron	Isopropyl alcohol 100%	14134722	1
Caelo	Isopropyl alcohol 100%	14134720	2
Caelo	Isopropyl alcohol 100%	14134724	1
Sanacorp	Isopropyl alcohol 100%	14134725	1
Caelo	Isopropyl alcohol 100%	14219204	2
Caelo	Isopropyl alcohol 100%	14219202	1
Caelo	Isopropyl alcohol 100%	14219205	1
Caelo	Isopropyl alcohol 100%	14219209	2
Caelo	Isopropyl alcohol 100%	14219206	1
Phoenix	Isopropyl alcohol 100%	14219204	1
Caelo	Isopropyl alcohol 100%	14219211	1
Caelo	Isopropyl alcohol 100%	14219213	2
Caelo	Isopropyl alcohol 100%	14219215	2
Caelo	Isopropyl alcohol 100%	14219214	2
Caelo	Isopropyl alcohol 100%	14219218	3
Caelo	Isopropyl alcohol 100%	14219220	1
Caelo	Isopropyl alcohol 100%	14219222	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	14219224	1
Caelo	Isopropyl alcohol 100%	14219228	2
Caelo	Isopropyl alcohol 100%	14219229	1
Caelo	Isopropyl alcohol 100%	14278901	3
noweda	Isopropyl alcohol 100%	14278905	1
Caelo	Isopropyl alcohol 100%	14278907	1
Caelo	Isopropyl alcohol 100%	14278911	2
Caelo	Isopropyl alcohol 100%	14278909	1
Sanacorp	Isopropyl alcohol 100%	14278913	1
Caelo	Isopropyl alcohol 100%	14278925	2
Caelo	Isopropyl alcohol 100%	14278921	1
Fagron	Isopropyl alcohol 100%	14278927	1
Caelo	Isopropyl alcohol 100%	14278927	2
Caelo	Isopropyl alcohol 100%	14278929	1
Caelo	Isopropyl alcohol 100%	14369301	4
Caelo	Isopropyl alcohol 100%	14369306	1
Caelo	Isopropyl alcohol 100%	14369310	1
Caelo	Isopropyl alcohol 100%	14369307	3
Caelo	Isopropyl alcohol 100%	14369303	1
Sanacorp	Isopropyl alcohol 100%	14369316	2
Caelo	Isopropyl alcohol 100%	14369317	2
Caelo	Isopropyl alcohol 100%	14369322	1
Caelo	Isopropyl alcohol 100%	14369319	1
Caelo	Isopropyl alcohol 100%	14369329	1
Caelo	Isopropyl alcohol 100%	14369326	1
Caelo	Isopropyl alcohol 100%	14369333	1
Caelo	Isopropyl alcohol 100%	14369337	1
Caelo	Isopropyl alcohol 100%	14369334	1
Caelo	Isopropyl alcohol 100%	15057505	2
Caelo	Isopropyl alcohol 100%	15057503	1
Caelo	Isopropyl alcohol 100%	15057507	2
Caelo	Isopropyl alcohol 100%	15057510	1
Caelo	Isopropyl alcohol 100%	15057511	2
Caelo	Isopropyl alcohol 100%	15057508	1
Caelo	Isopropyl alcohol 100%	15057513	1
Caelo	Isopropyl alcohol 100%	15057520	1
Caelo	Isopropyl alcohol 100%	15057524	1
Caelo	Isopropyl alcohol 100%	15057528	2
Caelo	Isopropyl alcohol 100%	15057526	1
Caelo	Isopropyl alcohol 100%	15057518	1
Caelo	Isopropyl alcohol 100%	15057514	1
Caelo	Isopropyl alcohol 100%	15057515	1
Caelo	Isopropyl alcohol 100%	15111401	1
Hedinger/GEHE	Isopropyl alcohol 100%	1514Q1079	1
Caelo	Isopropyl alcohol 100%	15171301	1
Sanacorp	Isopropyl alcohol 100%	15171302	1
Caelo	Isopropyl alcohol 100%	15171303	6
Caelo	Isopropyl alcohol 100%	15171305	1
Caelo	Isopropyl alcohol 100%	15171306	1
Caelo	Isopropyl alcohol 100%	15171307	1
Caelo	Isopropyl alcohol 100%	15171311	2
Caelo	Isopropyl alcohol 100%	15171312	2
Caelo	Isopropyl alcohol 100%	15171313	1
Caelo	Isopropyl alcohol 100%	15171315	1
Caelo	Isopropyl alcohol 100%	15171314	2
Caelo	Isopropyl alcohol 100%	15171328	3
Caelo	Isopropyl alcohol 100%	15171327	1
Caelo	Isopropyl alcohol 100%	15171332	1
Caelo	Isopropyl alcohol 100%	15171331	1
Caelo	Isopropyl alcohol 100%	15171333	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	15171340	1
Caelo	Isopropyl alcohol 100%	15171341	3
Caelo	Isopropyl alcohol 100%	15171335	1
Caelo	Isopropyl alcohol 100%	15171342	2
Caelo	Isopropyl alcohol 100%	1517308	1
Sanacorp	Isopropyl alcohol 100%	15300902	1
Caelo	Isopropyl alcohol 100%	15300903	2
Caelo	Isopropyl alcohol 100%	15300904	2
Caelo	Isopropyl alcohol 100%	15300913	1
Caelo	Isopropyl alcohol 100%	15300912	1
Caelo	Isopropyl alcohol 100%	15300917	1
Caelo	Isopropyl alcohol 100%	15300918	1
Caelo	Isopropyl alcohol 100%	15300915	1
Caelo	Isopropyl alcohol 100%	15300921	1
Caelo	Isopropyl alcohol 100%	15300927	2
Caelo	Isopropyl alcohol 100%	15300926	1
Caelo	Isopropyl alcohol 100%	15300929	1
Caelo	Isopropyl alcohol 100%	15300931	1
Caelo	Isopropyl alcohol 100%	15398803	2
Caelo	Isopropyl alcohol 100%	15398804	1
Caelo	Isopropyl alcohol 100%	15398806	3
Caesar&Loretz/Noweda	Isopropyl alcohol 100%	15398808	1
Caelo	Isopropyl alcohol 100%	15398808	1
Caelo	Isopropyl alcohol 100%	15398814	1
Caelo	Isopropyl alcohol 100%	15398811	1
Caelo	Isopropyl alcohol 100%	15398816	3
Caelo	Isopropyl alcohol 100%	15398819	1
Caelo	Isopropyl alcohol 100%	15398826	1
Caelo	Isopropyl alcohol 100%	15398820	1
Caelo	Isopropyl alcohol 100%	15398831	2
Caelo	Isopropyl alcohol 100%	15498831	1
Caelo	Isopropyl alcohol 100%	16039901	2
Sanacorp	Isopropyl alcohol 100%	16039901	1
Caelo	Isopropyl alcohol 100%	16039903	1
Caelo	Isopropyl alcohol 100%	16039905	1
Caelo	Isopropyl alcohol 100%	16039902	1
Caelo	Isopropyl alcohol 100%	16039904	2
Caelo	Isopropyl alcohol 100%	16039906	1
Caelo	Isopropyl alcohol 100%	16039907	1
Caelo	Isopropyl alcohol 100%	16039913	1
Caelo	Isopropyl alcohol 100%	16039912	1
Caelo	Isopropyl alcohol 100%	16039915	1
Caelo	Isopropyl alcohol 100%	16039921	1
Caelo	Isopropyl alcohol 100%	16039922	2
Caelo	Isopropyl alcohol 100%	16039926	1
Caelo	Isopropyl alcohol 100%	16039930	1
Caelo	Isopropyl alcohol 100%	16039925	1
Caelo	Isopropyl alcohol 100%	16039932	1
Caelo	Isopropyl alcohol 100%	1603997	1
Caelo	Isopropyl alcohol 100%	16137301	3
Noweda	Isopropyl alcohol 100%	16137302	1
Caelo	Isopropyl alcohol 100%	16137304	2
Caelo	Isopropyl alcohol 100%	16137303	1
Caelo	Isopropyl alcohol 100%	16137305	1
Caelo	Isopropyl alcohol 100%	16137322	2
Caelo	Isopropyl alcohol 100%	16137312	2
Caelo	Isopropyl alcohol 100%	16137311	1
Caelo	Isopropyl alcohol 100%	16137308	2
noweda	Isopropyl alcohol 100%	16137324	1
Caelo	Isopropyl alcohol 100%	16137327	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	16137325/161123AI	1
Caelo	Isopropyl alcohol 100%	16137332	1
Caelo	Isopropyl alcohol 100%	16137330	1
Caelo	Isopropyl alcohol 100%	16137336	3
Caelo	Isopropyl alcohol 100%	16137334	1
Caelo	Isopropyl alcohol 100%	16234001	2
Caelo	Isopropyl alcohol 100%	16234002	1
Caelo	Isopropyl alcohol 100%	16234004	2
Caelo	Isopropyl alcohol 100%	16234005	2
Caelo	Isopropyl alcohol 100%	16234010	3
Caelo	Isopropyl alcohol 100%	16234011	1
Caelo	Isopropyl alcohol 100%	16234014	1
Caelo	Isopropyl alcohol 100%	16234017	2
Caelo	Isopropyl alcohol 100%	16234018	1
Caelo	Isopropyl alcohol 100%	16234019	1
Caelo	Isopropyl alcohol 100%	16234025	2
Caelo	Isopropyl alcohol 100%	16234026	1
Caelo	Isopropyl alcohol 100%	16310601	1
Caesar /phönix	Isopropyl alcohol 100%	16310603	1
Caelo	Isopropyl alcohol 100%	16310605	2
Caelo	Isopropyl alcohol 100%	16310607	2
Caelo	Isopropyl alcohol 100%	16234031	1
Caelo	Isopropyl alcohol 100%	16310609	1
Caelo	Isopropyl alcohol 100%	16310611	1
Caelo	Isopropyl alcohol 100%	16310614	1
Hedinger GmbH	Isopropyl alcohol 100%	16310612	1
Caelo	Isopropyl alcohol 100%	16310612	1
Caelo	Isopropyl alcohol 100%	16310621	1
Caelo	Isopropyl alcohol 100%	16310623	1
Caelo	Isopropyl alcohol 100%	16310633	2
Henry Lamotte Oils	Isopropyl alcohol 100%	16310633	1
Caelo	Isopropyl alcohol 100%	16310634	1
Caelo	Isopropyl alcohol 100%	17044001	1
Caelo	Isopropyl alcohol 100%	17044006	1
Caelo	Isopropyl alcohol 100%	163234002	1
Caelo	Isopropyl alcohol 100%	17044003	1
Caelo	Isopropyl alcohol 100%	17044008	1
Caelo	Isopropyl alcohol 100%	17044002	1
Caesar & Loretz/Noweda	Isopropyl alcohol 100%	17044009	1
Caelo	Isopropyl alcohol 100%	17044013	1
Caelo	Isopropyl alcohol 100%	17044022	1
Caelo	Isopropyl alcohol 100%	17044028	1
Caelo	Isopropyl alcohol 100%	17189308	1
Caelo	Isopropyl alcohol 100%	17189313	2
Caelo	Isopropyl alcohol 100%	17189317	1
Caelo	Isopropyl alcohol 100%	17189319	2
Caelo	Isopropyl alcohol 100%	17300415	1
Caelo	Isopropyl alcohol 100%	17189321	1
Caelo	Isopropyl alcohol 100%	17300433	1
Caelo	Isopropyl alcohol 100%	17300432	1
Caelo	Isopropyl alcohol 100%	180109AI/17189332	1
Caelo	Isopropyl alcohol 100%	19051503	1
Caelo	Isopropyl alcohol 100%	19051405	1
Caelo	Isopropyl alcohol 100%	20011703	1
Caelo	Isopropyl alcohol 100%	21111402	1
Caelo	Isopropyl alcohol 100%	25021501	1
Caelo	Isopropyl alcohol 100%	27031502	1
Caelo	Isopropyl alcohol 100%	27091401	1
Caelo	Isopropyl alcohol 100%	29071401	1
Caelo	Isopropyl alcohol 100%	29011505	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	30041508	1
sophien-arzneimittel	Isopropyl alcohol 100%	310190514069	1
Caelo	Isopropyl alcohol 100%	31081501	1
Fischar	Isopropyl alcohol 100%	4105104	1
Fischar/Phönix	Isopropyl alcohol 100%	4103112	1
Hedinger/Noweda	Isopropyl alcohol 100%	4716M-01079	1
Hedinger/Noweda	Isopropyl alcohol 100%	4916U-01079	1
Hedinger	Isopropyl alcohol 100%	70587B	2
Caelo	Isopropyl alcohol 100%	6324064	1
Phönix 07.01.15	Isopropyl alcohol 100%	70650B	1
Sanacorp	Isopropyl alcohol 100%		1
hedinger	Isopropyl alcohol 100%	105208	1
Caelo	Isopropyl alcohol 100%	2011403	1
Caelo	Isopropyl alcohol 100%	1071301	1
Hedinger	Isopropyl alcohol 100%	231307	1
Hedinger	Isopropyl alcohol 100%	402211079	1
Hedinger	Isopropyl alcohol 100%	406211079	1
Caelo	Isopropyl alcohol 100%	4091201	1
Caelo	Isopropyl alcohol 100%	6061320	1
Caelo	Isopropyl alcohol 100%	5031303	1
Hedinger	Isopropyl alcohol 100%	50700B	1
Caelo	Isopropyl alcohol 100%	6071201	1
Caelo	Isopropyl alcohol 100%	6081204	1
Caelo	Isopropyl alcohol 100%	6081303	1
Hersteller: Hedinger	Isopropyl alcohol 100%	70236B	1
Hedinger	Isopropyl alcohol 100%	70416B	1
Kehr	Isopropyl alcohol 100%	70416B	1
Gehe	Isopropyl alcohol 100%	70382B	1
Hedinger	Isopropyl alcohol 100%	70219B	1
Gehe	Isopropyl alcohol 100%	704098	1
Hedinger	Isopropyl alcohol 100%	70422B	1
Hedinger/Anzag 27.02.2012	Isopropyl alcohol 100%	70417B	1
HedingerStuttgart/GEHE	Isopropyl alcohol 100%	70430B	1
Caelo	Isopropyl alcohol 100%	70439B	1
Hedinger	Isopropyl alcohol 100%	70439B	1
Hedinger über Noweda	Isopropyl alcohol 100%	70439B	1
Hedinger	Isopropyl alcohol 100%	70447B	2
Kehr	Isopropyl alcohol 100%	70439B	1
Hedinger/ Gehe	Isopropyl alcohol 100%	70444B	1
Hedinger	Isopropyl alcohol 100%	70448B	1
Hedinger/Holdermann 30...	Isopropyl alcohol 100%	70444B	1
Hedinger	Isopropyl alcohol 100%	70444B	1
Hedinger / Apotheke E.B...	Isopropyl alcohol 100%	70448B	1
Hedinger	Isopropyl alcohol 100%	70459B	2
Noweda	Isopropyl alcohol 100%	70459B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70457B	1
Hedinger	Isopropyl alcohol 100%	70461B	3
Hedinger	Isopropyl alcohol 100%	70457B	1
Hedinger / Apotheke E.B...	Isopropyl alcohol 100%	70461B	1
Hedinger	Isopropyl alcohol 100%	70463B	2
Hedinger	Isopropyl alcohol 100%	704638	1
Hedinger/Noweda	Isopropyl alcohol 100%	70465B	1
Hedinger	Isopropyl alcohol 100%	70465B	1
Hedinger, Heiligenwiese...	Isopropyl alcohol 100%	70471B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70471B	2
Hedinger/Anzag	Isopropyl alcohol 100%	70471B	1
Hedinger	Isopropyl alcohol 100%	70471B	2
Noweda	Isopropyl alcohol 100%	70475B	2
Fiebig 09.11.12 EK28,33	Isopropyl alcohol 100%	704758	1
Hedinger/Anzag	Isopropyl alcohol 100%	70475B	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Dr.Seeger GmbH	Isopropyl alcohol 100%	70475B	1
Klenk	Isopropyl alcohol 100%	70475b	1
Noweda/Hedinger	Isopropyl alcohol 100%	70476B	2
Jenne, Hedinger	Isopropyl alcohol 100%	70475B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70476B	2
Hedinger	Isopropyl alcohol 100%	70478B	5
Hedinger/VDL	Isopropyl alcohol 100%	704788	1
Hedinger/Gehe	Isopropyl alcohol 100%	70478B	1
Caelo	Isopropyl alcohol 100%	70478B	1
Hedinger Jenne	Isopropyl alcohol 100%	70478B	1
Kehr	Isopropyl alcohol 100%	70478B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70478B	3
Hedinger	Isopropyl alcohol 100%	70481B	2
Hedinger/Gehe	Isopropyl alcohol 100%	70481B	1
sanacorp	Isopropyl alcohol 100%	70481B	1
Hedinger	Isopropyl alcohol 100%	704888	1
Aug.Hedinger/Noweda	Isopropyl alcohol 100%	70488B	1
Hedinger	Isopropyl alcohol 100%	70488B	2
Hedinger /Sana	Isopropyl alcohol 100%	70488B	1
Noweda	Isopropyl alcohol 100%	70490B	1
Caelo	Isopropyl alcohol 100%	70492B	2
Hedinger/Gehe	Isopropyl alcohol 100%	70492B	1
Hedinger	Isopropyl alcohol 100%	70495B	1
HedingerStuttgart/GEHE	Isopropyl alcohol 100%	704958	1
Gehe 24.05.2013	Isopropyl alcohol 100%	70495B	1
HedingerStuttgart/GEHE	Isopropyl alcohol 100%	704978	1
Hedinger	Isopropyl alcohol 100%	70498B	1
Hedinger	Isopropyl alcohol 100%	70497B	1
Caelo	Isopropyl alcohol 100%	70498B	3
Anzag 08.02.2013	Isopropyl alcohol 100%	70500B	1
Hedinger	Isopropyl alcohol 100%	70500B	1
Caesar &Loretz GmbH	Isopropyl alcohol 100%	70503B	1
Hedinger	Isopropyl alcohol 100%	70501B	4
Hedinger/Gehe	Isopropyl alcohol 100%	70501B	1
Hedinger/Alliance Healt...	Isopropyl alcohol 100%	70503B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70505B	1
Hedinger	Isopropyl alcohol 100%	70507B	7
Hedinger/Gehe	Isopropyl alcohol 100%	70507B	1
Gehe	Isopropyl alcohol 100%	70507B	1
Hedinger/AHD	Isopropyl alcohol 100%	70507B	1
Hedinger	Isopropyl alcohol 100%	70511B	3
Hedinger/ Noweda	Isopropyl alcohol 100%	70517B	1
Bombastus	Isopropyl alcohol 100%	70517B	1
Hedinger, Heiligenwiese...	Isopropyl alcohol 100%	70519B	1
Dr.Seeger GmbH	Isopropyl alcohol 100%	70519B	1
Hedinger/Anzag	Isopropyl alcohol 100%	705198	1
Hedinger/Noweda	Isopropyl alcohol 100%	70519B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70519B	1
Noweda	Isopropyl alcohol 100%	70519B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70520B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70519B	1
Hedinger/ Noweda	Isopropyl alcohol 100%	70520B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70525B	1
Hedinger/Jenne	Isopropyl alcohol 100%	70525B	1
Hedinger	Isopropyl alcohol 100%	70527B	6
Stumpf/Hedinger	Isopropyl alcohol 100%	70527B	2
Dr.Seeger GmbH	Isopropyl alcohol 100%	70527B	1
HedingerStuttgart/GEHE	Isopropyl alcohol 100%	70527B	1
Anzag/Seeger	Isopropyl alcohol 100%	70527B	1
Klenk	Isopropyl alcohol 100%	70527B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Isopropyl alcohol 100%	70530B	2
Alliance/ Hedinger	Isopropyl alcohol 100%	70530B	3
Hedinger/Gehe	Isopropyl alcohol 100%	70527B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70532B	1
Hedinger	Isopropyl alcohol 100%	70532B	3
Hedinger/AHD	Isopropyl alcohol 100%	70532B	1
Phönix /Hedinger	Isopropyl alcohol 100%	70532B	2
Caelo	Isopropyl alcohol 100%	70536B	2
Hedinger	Isopropyl alcohol 100%	70536B	2
Phönix	Isopropyl alcohol 100%	70536B	1
Euro OTC	Isopropyl alcohol 100%	70536B	1
Hedinger Phönix	Isopropyl alcohol 100%	70536B	1
Alliance Healthcare/Hed...	Isopropyl alcohol 100%	70536B	1
Caelo	Isopropyl alcohol 100%	70540B	1
Stumpf	Isopropyl alcohol 100%	70540B	1
Hedinger/Gehe	Isopropyl alcohol 100%	70540B	1
Hedinger	Isopropyl alcohol 100%	70540B	1
Gehe/Hedinger	Isopropyl alcohol 100%	70540B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70540B	1
Alliance/ Hedinger	Isopropyl alcohol 100%	70540B	1
Fagron	Isopropyl alcohol 100%	70541	1
Noweda/ Dr. Seeger	Isopropyl alcohol 100%	70540B	1
Caelo	Isopropyl alcohol 100%	70541B	2
Gehe	Isopropyl alcohol 100%	70541B	1
Hedinger / Apotheke E.B...	Isopropyl alcohol 100%	70541B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70541B	1
Anzag/Hedinger	Isopropyl alcohol 100%	70541B	6
Hedinger	Isopropyl alcohol 100%	70541B	4
Hedinger/ jenne	Isopropyl alcohol 100%	70541B	1
Hedinger	Isopropyl alcohol 100%	70545B	4
Caelo	Isopropyl alcohol 100%	70545B	1
Hedinger/Anzag	Isopropyl alcohol 100%	70545B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70545B	1
Hedinger/Phoenix	Isopropyl alcohol 100%	70545B	2
Apoth.E.Bartel	Isopropyl alcohol 100%	70545B	2
Hedinger/ Sanacorp	Isopropyl alcohol 100%	70545B	1
Hedinger/Noweda	Isopropyl alcohol 100%	7054907	1
Hedinger	Isopropyl alcohol 100%	70545b	1
Dr.Seeger GmbH	Isopropyl alcohol 100%	70553B	2
Caelo	Isopropyl alcohol 100%	70553B	1
Alliance/ Hedinger	Isopropyl alcohol 100%	70553B	2
Hedinger/Gehe	Isopropyl alcohol 100%	70553B	1
Fischar	Isopropyl alcohol 100%	70553B	1
Hedinger	Isopropyl alcohol 100%	70553B	5
Hedinger/Anzag	Isopropyl alcohol 100%	70553B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70553BH	1
Hedinger	Isopropyl alcohol 100%	705540B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70555B	2
Caelo	Isopropyl alcohol 100%	70555B	1
Hedinger	Isopropyl alcohol 100%	70555B	4
Hedinger/ Noweda	Isopropyl alcohol 100%	70555B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70555B	1
Hedinger/Jenne	Isopropyl alcohol 100%	70559B	1
Alliance/ Hedinger	Isopropyl alcohol 100%	70559B	1
Hedinger	Isopropyl alcohol 100%	70559B	1
Hedinger/Sanacorp	Isopropyl alcohol 100%	70564B	1
Anzag/Hedinger	Isopropyl alcohol 100%	70562B	1
Caelo	Isopropyl alcohol 100%	70564B	1
Hedinger , Stuttgart	Isopropyl alcohol 100%	70565B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70565B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Noweda	Isopropyl alcohol 100%	70565B	1
Caelo	Isopropyl alcohol 100%	70565B	1
Hedinger	Isopropyl alcohol 100%	705728	2
Noweda/Hedinger	Isopropyl alcohol 100%	70572B	1
Hedinger/Anzag	Isopropyl alcohol 100%	705728	1
Hedinger	Isopropyl alcohol 100%	705738	2
Hedinger	Isopropyl alcohol 100%	70574B	3
Hedinger	Isopropyl alcohol 100%	70573B	1
HEDINGER / 20	Isopropyl alcohol 100%	70573B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70574B	1
Alliance Healthcare/Hed...	Isopropyl alcohol 100%	70577B	1
Hedinger	Isopropyl alcohol 100%	70577B	1
Hedinger	Isopropyl alcohol 100%	70578	1
Anzag	Isopropyl alcohol 100%	70578B	2
Hedinger/Sanacorp	Isopropyl alcohol 100%	70578B	2
Hedinger	Isopropyl alcohol 100%	70578B	2
Euro OTC	Isopropyl alcohol 100%	70578B	1
Hedinger/Otto	Isopropyl alcohol 100%	70581B	1
Hedinger/Noweda	Isopropyl alcohol 100%	70584B	1
Hedinger	Isopropyl alcohol 100%	70581B	1
Sanacorp	Isopropyl alcohol 100%	70593B	1
Caelo	Isopropyl alcohol 100%	70593B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70584B	1
Hedinger / Noweda	Isopropyl alcohol 100%	70587B	1
hedinger	Isopropyl alcohol 100%	70584B	1
Hedinger	Isopropyl alcohol 100%	70595B	1
Anzag	Isopropyl alcohol 100%	7059B	1
Noweda/Hedinger	Isopropyl alcohol 100%	70610B	1
Noweda/Hedinger	Isopropyl alcohol 100%	706148	1
Hedinger	Isopropyl alcohol 100%	70614B	1
Caelo	Isopropyl alcohol 100%	70701B	1
Hedinger / Sanacorp	Isopropyl alcohol 100%	707446B	1
Caelo	Isopropyl alcohol 100%	713Q-01079	1
Caelo	Isopropyl alcohol 100%	9101314	1
Caelo	Isopropyl alcohol 100%	8012014E	1
Caelo	Isopropyl alcohol 100%	750520B	1
Caelo	Isopropyl alcohol 100%	10121304	1
Caelo	Isopropyl alcohol 100%	10051302	1
Caelo	Isopropyl alcohol 100%	11371204	1
Caelo	Isopropyl alcohol 100%	108073	1
Caelo	Isopropyl alcohol 100%	11311602	1
Noweda	Isopropyl alcohol 100%	11371222	1
Caelo	Isopropyl alcohol 100%	11371214	1
Caelo	Isopropyl alcohol 100%	11371212	1
Caelo	Isopropyl alcohol 100%	11371217	1
Caelo	Isopropyl alcohol 100%	12031106	1
Caelo	Isopropyl alcohol 100%	12031129	1
Gehe	Isopropyl alcohol 100%	12031105	1
Hedinger	Isopropyl alcohol 100%	1207211079	1
Caelo	Isopropyl alcohol 100%	12091305	1
Caelo	Isopropyl alcohol 100%	12112206	1
Caelo	Isopropyl alcohol 100%	12112220	1
Caelo	Isopropyl alcohol 100%	12112224	3
Noweda	Isopropyl alcohol 100%	12112218	1
Caelo	Isopropyl alcohol 100%	12112209	1
Noweda	Isopropyl alcohol 100%	12112213	1
Caelo	Isopropyl alcohol 100%	12179307	3
Großhandel	Isopropyl alcohol 100%	12179309	1
Caelo	Isopropyl alcohol 100%	12179310	1
Caelo	Isopropyl alcohol 100%	12179311	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	12179321	1
Noweda	Isopropyl alcohol 100%	12179322	1
Caelo	Isopropyl alcohol 100%	12179317	2
Caelo	Isopropyl alcohol 100%	12179326	1
Caelo	Isopropyl alcohol 100%	12254702	1
Caelo	Isopropyl alcohol 100%	12179328	1
Caelo	Isopropyl alcohol 100%	12254708	1
Caelo	Isopropyl alcohol 100%	12254706	1
Caelo	Isopropyl alcohol 100%	12254708/129	1
Caelo	Isopropyl alcohol 100%	12254709	2
Gehe	Isopropyl alcohol 100%	12254718	1
Caelo	Isopropyl alcohol 100%	12254714	1
Noweda	Isopropyl alcohol 100%	12254715	1
Caelo	Isopropyl alcohol 100%	12254722	1
Großhandel	Isopropyl alcohol 100%	12254719	1
Noweda	Isopropyl alcohol 100%	12254724	1
Caelo	Isopropyl alcohol 100%	12254719	1
Caelo	Isopropyl alcohol 100%	12254724	1
Caelo	Isopropyl alcohol 100%	12254728	1
Caelo	Isopropyl alcohol 100%	12254729	1
Caelo	Isopropyl alcohol 100%	12321004	1
Phönix 11.12.12	Isopropyl alcohol 100%	12321005	1
Caelo	Isopropyl alcohol 100%	12321008	1
Caelo	Isopropyl alcohol 100%	12321012	1
Caelo	Isopropyl alcohol 100%	12321015	1
Caelo	Isopropyl alcohol 100%	12321016	1
Caelo	Isopropyl alcohol 100%	12321022	1
Caelo	Isopropyl alcohol 100%	12321020	2
Caelo	Isopropyl alcohol 100%	12321023	3
Phönix/ henry lamotte	Isopropyl alcohol 100%	12321024	1
Caelo	Isopropyl alcohol 100%	12321031	1
Großhandel	Isopropyl alcohol 100%	12321026	1
Caelo	Isopropyl alcohol 100%	12321032	1
Großhandel	Isopropyl alcohol 100%	12321029	1
Anzag	Isopropyl alcohol 100%	12321028	1
Caelo	Isopropyl alcohol 100%	12379501	1
Caelo	Isopropyl alcohol 100%	12379502	3
Caelo	Isopropyl alcohol 100%	12379508	2
Caelo	Isopropyl alcohol 100%	12379510	1
Caelo	Isopropyl alcohol 100%	12379513	1
Noweda	Isopropyl alcohol 100%	12379510	1
Noweda	Isopropyl alcohol 100%	12379515	1
Caelo	Isopropyl alcohol 100%	12379515	2
Caelo	Isopropyl alcohol 100%	12379517	1
Caelo	Isopropyl alcohol 100%	12379518	2
Caelo	Isopropyl alcohol 100%	12379520	2
Caelo	Isopropyl alcohol 100%	12379521	2
Caelo	Isopropyl alcohol 100%	12379522	1
Caelo	Isopropyl alcohol 100%	12379525	2
Großhandel	Isopropyl alcohol 100%	12379525	1
Phönix 15.07.2013	Isopropyl alcohol 100%	12379528	1
Caelo	Isopropyl alcohol 100%	12379527	1
Caelo	Isopropyl alcohol 100%	13021409	1
Caelo	Isopropyl alcohol 100%	13050904	2
Caelo	Isopropyl alcohol 100%	13050902	1
Caelo	Isopropyl alcohol 100%	13050901	1
Caelo	Isopropyl alcohol 100%	13050907	2
Erika apotheke	Isopropyl alcohol 100%	13050911	1
Caelo	Isopropyl alcohol 100%	13050911	1
Caelo	Isopropyl alcohol 100%	13050914	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 100%	13050916	3
Caelo	Isopropyl alcohol 100%	13050917	2
Caelo	Isopropyl alcohol 100%	13050920	1
Caelo	Isopropyl alcohol 100%	13050925	2
Caelo	Isopropyl alcohol 100%	13050922	1
Caelo	Isopropyl alcohol 100%	13050928	1
Caelo	Isopropyl alcohol 100%	13050927	1
Caelo	Isopropyl alcohol 100%	13071304	1
Caelo	Isopropyl alcohol 100%	13050929	2
Caelo	Isopropyl alcohol 100%	13148003	1
Caelo	Isopropyl alcohol 100%	13148012	1
Caelo	Isopropyl alcohol 100%	13148010	1
Caelo	Isopropyl alcohol 100%	13148015	3
Caelo	Isopropyl alcohol 100%	13148005	1
Caelo	Isopropyl alcohol 100%	13148016	1
Caelo	Isopropyl alcohol 100%	13148019	3
Caelo	Isopropyl alcohol 100%	13148024	1
Großhandel	Isopropyl alcohol 100%	13148027	1
Caelo	Isopropyl alcohol 100%	13148021	1
Caelo	Isopropyl alcohol 100%	13148028	1
Sanacorp	Isopropyl alcohol 100%	1324919	2
Caelo	Isopropyl alcohol 100%	132/AKD54	1
Caelo	Isopropyl alcohol 100%	1324931	1
Caelo	Isopropyl alcohol 100%	13249403	2
Caelo	Isopropyl alcohol 100%	13249401	1
Gehe	Isopropyl alcohol 100%	13249403	1
Euro OTC	Isopropyl alcohol 100%	13249405	1
Caelo	Isopropyl alcohol 100%	13249404	1
Caelo	Isopropyl alcohol 100%	13249405	1
Caelo	Isopropyl alcohol 100%	13249407	1
Caelo	Isopropyl alcohol 100%	13249412	2
Caelo	Isopropyl alcohol 100%	13249411	1
Caelo	Isopropyl alcohol 100%	13249413	2
Caelo	Isopropyl alcohol 100%	13249414	2
Caelo	Isopropyl alcohol 100%	13249415	3
Caelo	Isopropyl alcohol 100%	13249424	1
Caelo	Isopropyl alcohol 100%	13249416	1
Sanacorp	Isopropyl alcohol 100%	13249416	1
Caelo	Isopropyl alcohol 100%	13249427	1
Caelo	Isopropyl alcohol 100%	13249430	1
Caelo	Isopropyl alcohol 100%	13249423	1
Caelo	Isopropyl alcohol 100%	13249420	1
Caelo	Isopropyl alcohol 100%	13249429	1
Phönix	Isopropyl alcohol 100%	13249432	1
Caelo	Isopropyl alcohol 100%	13249435	3
Caelo	Isopropyl alcohol 100%	13249434	1
Caelo	Isopropyl alcohol 100%	13336911	1
Caelo	Isopropyl alcohol 100%	13336912	2
Caelo	Isopropyl alcohol 100%	13249437	1
Caelo	Isopropyl alcohol 100%	1329423	1
Caelo	Isopropyl alcohol 100%	13336915	1
Caelo	Isopropyl alcohol 100%	1329421	1
Caelo	Isopropyl alcohol 100%	13336930	1
Caelo	Isopropyl alcohol 100%	13336924	1
Caelo	Isopropyl alcohol 100%	13336922	1
Caelo	Isopropyl alcohol 100%	13336931	1
Caelo	Isopropyl alcohol 100%	13336933	1
Caelo	Isopropyl alcohol 100%	13336934	1
Caelo	Isopropyl alcohol 100%	13336938	3
Großhandel	Isopropyl alcohol 100%	13336935	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
cealo	Isopropyl alcohol 100%	13336936	1
Caelo	Isopropyl alcohol 100%	13336936	1
Caelo	Isopropyl alcohol 100%	13336939	2
Caelo	Isopropyl alcohol 100%	13336941	2
Caelo	Isopropyl alcohol 100%	13336942	1
Caelo	Isopropyl alcohol 100%	13336944	1
Caelo	Isopropyl alcohol 100%	13336943	2
Caelo	Isopropyl alcohol 100%	13397304	2
Caelo	Isopropyl alcohol 100%	13336945	1
Noweda	Isopropyl alcohol 100%	13397303	1
Caelo	Isopropyl alcohol 100%	13397311	1
Caelo	Isopropyl alcohol 100%	13397323	1
Caelo	Isopropyl alcohol 100%	140510265	1
Fagron	Isopropyl alcohol 100%	14134702	1
Phonix	Isopropyl alcohol 100%	14134715	1
Caelo	Isopropyl alcohol 100%	14369370	1
Anzag	Isopropyl alcohol 100%	14716	1
Anzag	Isopropyl alcohol 100%	14219233	1
Hedinger	Isopropyl alcohol 100%	1512211079	1
Anzag	Isopropyl alcohol 100%	14717	1
Caelo	Isopropyl alcohol 100%	15300910	1
Caelo	Isopropyl alcohol 100%	15300930	1
Caelo	Isopropyl alcohol 100%	16011406	1
Caelo	Isopropyl alcohol 100%	15398813	1
Hedinger	Isopropyl alcohol 100%	1601311079	1
Caelo	Isopropyl alcohol 100%	16039919	1
Celo	Isopropyl alcohol 100%	16039925	1
Celo	Isopropyl alcohol 100%	16039937	1
Caelo	Isopropyl alcohol 100%	20081204	1
Hedinger	Isopropyl alcohol 100%	2x070461B	1
Caelo	Isopropyl alcohol 100%	18041306	1
Caelo	Isopropyl alcohol 100%	16137326	1
Caelo	Isopropyl alcohol 100%	21021311	1
Caelo	Isopropyl alcohol 100%	20060912-4	1
Caelo	Isopropyl alcohol 100%	2201201301	1
Caelo	Isopropyl alcohol 100%	22011302	1
Caelo	Isopropyl alcohol 100%	22091203	1
Caelo	Isopropyl alcohol 100%	22111204	1
Hedinger	Isopropyl alcohol 100%	2310211079	1
Caelo	Isopropyl alcohol 100%	24071201	1
Caelo	Isopropyl alcohol 100%	22121201	1
Caelo	Isopropyl alcohol 100%	25061201	1
Caelo	Isopropyl alcohol 100%	29101304	1
Caelo	Isopropyl alcohol 100%	29051202	1
Caelo	Isopropyl alcohol 100%	28021316	1
Caelo	Isopropyl alcohol 100%	30111303	1
Caelo	Isopropyl alcohol 100%	28081302	1
Hedinger	Isopropyl alcohol 100%	3013E1079	1
Kost/ - /16.09.13	Isopropyl alcohol 100%	310120913059	1
kost	Isopropyl alcohol 100%	310120913059	2
Caelo	Isopropyl alcohol 100%	31101204	1
Kost/kost/sophien	Isopropyl alcohol 100%	310261113062	1
Fischer/Phönix	Isopropyl alcohol 100%	4103014	1
Gehe 12.04.2013	Isopropyl alcohol 100%	4912Q-07001	1
Aug. Hedinger GmbH & Co...	Isopropyl alcohol 100%	501114B	1
Caelo	Isopropyl alcohol 100%	543892	1
Caelo	Isopropyl alcohol 100%	7272262	1
spinnenrad	Isopropyl alcohol 100%	611410093	1
hedinger/anzag	Isopropyl alcohol 100%	705838	1
Caelo	Isopropyl alcohol 100%	7411	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Isopropyl alcohol 100%		1
Hedinger/Otto	Isopropyl alcohol 100%		1
Hedinger	Isopropyl alcohol 100%	7474942	1
Caelo	Isopropyl alcohol 100%	7776	1

- 89 618 spectra from 501 *Apo-Ident* customers from a total of 26 242 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Isopropyl alcohol 100%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Isopropyl alcohol 100%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	140	0	93 734
Type C	0	1679	55	89 618

The substance/substance group *Isopropyl alcohol 100%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6586 %)	96.8281 % (> 96.6551 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30759	31517	5.89	203.48
31232	31517	1.47	197.88
31386	31517	3.04	195.73
31517	31517	0.00	196.53
32817	31517	6.05	186.07

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Isopropyl alcohol 70%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30848-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Isopropyl alcohol 70%; Alcohol isopropylus 70%

Special notes

When selecting the *Isopropyl alcohol 70%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Isopropyl alcohol 70%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Isopropyl alcoho...	14076503	31388	40	1412279
Caelo	Isopropyl alcoho...	14331004	31622	40	20141104
Caelo	Isopropyl alcoho...	15403802	31831	40	20151127
Caelo	Isopropyl alcoho...	163189	33105	40	20161206*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Isopropyl alcohol 70%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 760 spectra from 20 reference samples of the substance/substance group *Isopropyl alcohol 70%*.
- These include spectra of independent samples from 17 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
zu erfassen	Isopropyl alcohol 70%	zu erfassen	30391	40
Hofmann's	Isopropyl alcohol 70%	301012	30848	40
Caelo	Isopropyl alcohol 70%	14331004	31622 [†]	20
Caelo	Isopropyl alcohol 70%	15403802	31831 [†]	20
Caelo	Isopropyl alcohol 70%	160884	32216	40
Caelo	Isopropyl alcohol 70%	161678	32318	40
Caelo	Isopropyl alcohol 70%	161677	32319	40
Caelo	Isopropyl alcohol 70%	161676	32320	40
Caelo	Isopropyl alcohol 70%	161679	32322	40
Caelo	Isopropyl alcohol 70%	163190	33106	40
Caelo	Isopropyl alcohol 70%	163191	33107	40
Caelo	Isopropyl alcohol 70%	163192	33108	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Isopropyl alcohol 70%	172029	33685	40
Caelo	Isopropyl alcohol 70%	172028	33686	40
Caelo	Isopropyl alcohol 70%	172031	33689	40
Caelo	Isopropyl alcohol 70%	170821	33702	40
Caelo	Isopropyl alcohol 70%	181143	34067	40
Caelo	Isopropyl alcohol 70%	181142	34068	40
Caelo	Isopropyl alcohol 70%	181141	34069	40
Caelo	Isopropyl alcohol 70%	180739	34071	40

- 93114 spectra from a total of 2148 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1379 spectra from 180 *Apo-Ident* customers from 863 batches from the substance/substance group *Isopropyl alcohol 70%*.
- These include spectra of independent samples from 856 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Hedinger/Anz	Isopropyl alcohol 70%	53399	1
Alte Apo Defektur	Isopropyl alcohol 70%	117A-01080	1
Gehe	Isopropyl alcohol 70%	53391	1
Alte Apo Defektur	Isopropyl alcohol 70%	517I-01080	1
A. Hedinger/Noweda	Isopropyl alcohol 70%	53336	2
Hedinger/Noweda	Isopropyl alcohol 70%	53406	1
Hedinger/Phönix	Isopropyl alcohol 70%	53410	3
Hedinger/Noweda	Isopropyl alcohol 70%	53417	2
Sanacorp	Isopropyl alcohol 70%	53417	1
Hedinger	Isopropyl alcohol 70%	53425	1
Sanacorp	Isopropyl alcohol 70%	53430	1
GEHE	Isopropyl alcohol 70%	53420	2
Hedinger/VDL	Isopropyl alcohol 70%	53425	1
Hedinger/Fiebig	Isopropyl alcohol 70%	53440	1
Hedinger/Kehr	Isopropyl alcohol 70%	53435	1
Fiebig/Hedinger	Isopropyl alcohol 70%	53445	1
Hedinger/Fiebig	Isopropyl alcohol 70%	53456	3
Phoenix	Isopropyl alcohol 70%	53447	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53465	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53466	3
Hedinger/AHD	Isopropyl alcohol 70%	53479	1
GEHE	Isopropyl alcohol 70%	53491	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53494	2

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Fiebig	Isopropyl alcohol 70%	53495	1
Hedinger/Fiebig	Isopropyl alcohol 70%	53505	1
Hedinger/Fiebig	Isopropyl alcohol 70%	53510	1
Caelo	Isopropyl alcohol 70%	13363101	2
Caelo	Isopropyl alcohol 70%	15437601	7
Hustadt-Apotheke	Isopropyl alcohol 70%	160413/2	1
Caelo	Isopropyl alcohol 70%	16167902	1
Alte Apo Defektur	Isopropyl alcohol 70%	1617M-01080	1
Gehe	Isopropyl alcohol 70%	16318901	1
Caelo	Isopropyl alcohol 70%	16318901	1
Gehe	Isopropyl alcohol 70%	16319001	1
Caelo	Isopropyl alcohol 70%	16319201	2
Caelo	Isopropyl alcohol 70%	17202802	1
Caelo	Isopropyl alcohol 70%	17311403	1
Caelo	Isopropyl alcohol 70%	18073801	1
Apo Dunningen	Isopropyl alcohol 70%	1917A-01080	1
Eigenherstellung Hubert...	Isopropyl alcohol 70%	2014072221 (DF)	1
Fischar	Isopropyl alcohol 70%	2183	1
Alte Apo Defektur	Isopropyl alcohol 70%	2217Q-01080	1
Alte Apo Defektur	Isopropyl alcohol 70%	2217M-01080	1
Caelo	Isopropyl alcohol 70%	2314E-01080	1
Caelo	Isopropyl alcohol 70%	2415Q-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	3215M-07123	1
Caelo	Isopropyl alcohol 70%	3414E-01080	1
Wolf-Apotheke / Defektur	Isopropyl alcohol 70%	41017/50W1o	1
Apo Dunningen	Isopropyl alcohol 70%	3517A-01080	1
Fischar	Isopropyl alcohol 70%	4114064	1
Fischar / Anzag	Isopropyl alcohol 70%	4115073	1
Fagron	Isopropyl alcohol 70%	4114104	1
Otto Fischar GmbH & Co...	Isopropyl alcohol 70%	4115114	2
Phoenix	Isopropyl alcohol 70%	4116125053416	1
Fiebig	Isopropyl alcohol 70%	4116123	1
Phoenix	Isopropyl alcohol 70%	4116126	6
OttoFischar/Fiebig	Isopropyl alcohol 70%	41170	1
Otto Fischer/Fiebig	Isopropyl alcohol 70%	4117041	15
Otto Fischer/Fiebig	Isopropyl alcohol 70%	4117042	11
Otto Fischer/Phönix	Isopropyl alcohol 70%	4117043	20
Fischar/Fiebig	Isopropyl alcohol 70%	4117064	7
Fischar/Fiebig	Isopropyl alcohol 70%	4117065	12
Caelo	Isopropyl alcohol 70%	4118012	1
Caelo	Isopropyl alcohol 70%	4914A-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	4814E-07123	1
Eigen	Isopropyl alcohol 70%	5216E-01079	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-021116-2D	1
Caelo	Isopropyl alcohol 70%	Caelo	1
Winthrop	Isopropyl alcohol 70%	2600750T	1
Hedinger/Anzag	Isopropyl alcohol 70%	303008	1
Fischar/Sanacorp	Isopropyl alcohol 70%	4114044	1
Eigene Herstellung	Isopropyl alcohol 70%	AL310	1
Labor	Isopropyl alcohol 70%	7939	1
Wintrop AM GmbH	Isopropyl alcohol 70%	2600749T	1
Hubertus Apotheke	Isopropyl alcohol 70%	ChB120914	1
Eigenherstellung	Isopropyl alcohol 70%	9091115	1
Eigenherstellung	Isopropyl alcohol 70%	9261114	1
Eigenherstellung	Isopropyl alcohol 70%	9211215	1
Eigenherstellung	Isopropyl alcohol 70%	9270215	1
Eigenherstellung	Isopropyl alcohol 70%	9100414	1
Hedinger/AHD	Isopropyl alcohol 70%	53326	1
Hedinger/AHD	Isopropyl alcohol 70%	53407	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	140612014	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 70%	14229202	1
Caelo	Isopropyl alcohol 70%	15150202	4
Caelo	Isopropyl alcohol 70%	15192003	2
Caelo	Isopropyl alcohol 70%	15308201	1
Caelo	Isopropyl alcohol 70%	15404002	2
Caelo	Isopropyl alcohol 70%	16319202	1
Caelo	Isopropyl alcohol 70%	1915A-01080	1
Fischar	Isopropyl alcohol 70%	4114082	1
Eigenherstellung	Isopropyl alcohol 70%	9010813	1
Eigenherstellung	Isopropyl alcohol 70%	9230413	1
Hedinger	Isopropyl alcohol 70%	HA32	1
Eigenherstellung	Isopropyl alcohol 70%	9070313	1
Caelo	Isopropyl alcohol 70%	RH141030-0003	1
GEHE	Isopropyl alcohol 70%	221160/3	1
Sino Phyto	Isopropyl alcohol 70%	11370004	3
GEHE	Isopropyl alcohol 70%	12037601	1
Caelo	Isopropyl alcohol 70%	12136303	6
Eigenherstellung	Isopropyl alcohol 70%	1513A-01080	1
Defektur	Isopropyl alcohol 70%	21012013	1
Brugger	Isopropyl alcohol 70%	21112013	1
Eigenproduktion	Isopropyl alcohol 70%	2413E-01080	1
Eigenherstellung	Isopropyl alcohol 70%	3313Q-01080	1
GEHE	Isopropyl alcohol 70%	603117000	1
Wolf-Apotheke	Isopropyl alcohol 70%		5
N.N.	Isopropyl alcohol 70%		1
eigene Herstellung	Isopropyl alcohol 70%	91215	1
Hubertus Apotheke	Isopropyl alcohol 70%	103171	1
Hubertus Apotheke	Isopropyl alcohol 70%	103172	1
Hubertus Apotheke	Isopropyl alcohol 70%	103173	1
sittig	Isopropyl alcohol 70%	10815	1
Eigenherstellung	Isopropyl alcohol 70%	1-11122014	1
Apotheke	Isopropyl alcohol 70%	116E-07013	1
Caelo	Isopropyl alcohol 70%	2041501	1
Caelo	Isopropyl alcohol 70%	116M-01080	1
eigene Herstellung	Isopropyl alcohol 70%	1201-15	1
Hochhaus-Apotheke	Isopropyl alcohol 70%	117E-01080	1
Eigenherstellung	Isopropyl alcohol 70%	1-20082014	1
Eigenherstellung	Isopropyl alcohol 70%	1-23102014	1
Eigen/Defektur	Isopropyl alcohol 70%	216A-01080	1
Defektur	Isopropyl alcohol 70%	2/15	1
Eigen/Defektur	Isopropyl alcohol 70%	216M-01080	1
Eigenherstellung	Isopropyl alcohol 70%	216M-01080	1
Sonnen-Apotheke	Isopropyl alcohol 70%	3/16	1
Rondell Apotheke	Isopropyl alcohol 70%	245Q-07123	1
Hubertus Apotheke	Isopropyl alcohol 70%	31215	1
Eigenherstellung	Isopropyl alcohol 70%	316E	1
Caelo	Isopropyl alcohol 70%	3081605	1
Eigen Defektur	Isopropyl alcohol 70%	316Q-01080	1
Apotheke Dunningen	Isopropyl alcohol 70%	317A-01080	1
Eigenherstellung	Isopropyl alcohol 70%	4022016-01	1
Stifts-Apotheke	Isopropyl alcohol 70%	317E-01080	1
Hirsch Apotheke	Isopropyl alcohol 70%	415M-0108	1
Eigenherstellung	Isopropyl alcohol 70%	515E-01080	1
Eigenherstellung	Isopropyl alcohol 70%	517Q-01080	1
Caelo	Isopropyl alcohol 70%	5061501	1
Defekur	Isopropyl alcohol 70%	516Q-01080	1
eigene Herstellung	Isopropyl alcohol 70%	5210115	1
Hedinger/GEHE	Isopropyl alcohol 70%	53266	1
Gehe	Isopropyl alcohol 70%	53264	1
Hedinger / Phönix	Isopropyl alcohol 70%	53266	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
gehe	Isopropyl alcohol 70%	53270	1
Hedinger/AHD	Isopropyl alcohol 70%	53277	1
Alliance Health 17.03.2...	Isopropyl alcohol 70%	53279	2
Fagron	Isopropyl alcohol 70%	53279	1
Phoenix 17.03.2014	Isopropyl alcohol 70%	53280	1
Hedinger	Isopropyl alcohol 70%	53280	5
Hedinger	Isopropyl alcohol 70%	53284	1
Alliance Healthcare/Hed...	Isopropyl alcohol 70%	53285	3
GEHE am 17.06.2014	Isopropyl alcohol 70%	53286	1
Hedinger1	Isopropyl alcohol 70%	53290	1
Gehe	Isopropyl alcohol 70%	53290	1
Hedinger / Phönix	Isopropyl alcohol 70%	53290	1
Phönix/Hedinger	Isopropyl alcohol 70%	53295	2
AHD/Hedinger	Isopropyl alcohol 70%	53296	1
Hedinger	Isopropyl alcohol 70%	53291	1
Gehe	Isopropyl alcohol 70%	53296	1
Hedinger/AHD	Isopropyl alcohol 70%	53296	1
GEHE am 16.08.2014	Isopropyl alcohol 70%	53297	1
Hedinger / Noweda	Isopropyl alcohol 70%	53300	1
Hedinger	Isopropyl alcohol 70%	53297	1
Hedinger	Isopropyl alcohol 70%	53301	1
Hedinger/Gehe	Isopropyl alcohol 70%	53302	1
Hedinger	Isopropyl alcohol 70%	53303	1
Gehe	Isopropyl alcohol 70%	53303	1
Hedinger/Geilenkirchen	Isopropyl alcohol 70%	53306	2
Phönix 05.01.15	Isopropyl alcohol 70%	53303	1
Hedinger	Isopropyl alcohol 70%	53306	1
Hedinger	Isopropyl alcohol 70%	53308	1
Hedinger	Isopropyl alcohol 70%	53309	3
Hedinger / Phönix	Isopropyl alcohol 70%	53309	1
Gehe/Hedinger	Isopropyl alcohol 70%	53310	1
Hedinger / Noweda	Isopropyl alcohol 70%	5331	1
Hedinger/AHD	Isopropyl alcohol 70%	53310	1
Hedinger/Noweda	Isopropyl alcohol 70%	53310	1
Hedinger/Gehe	Isopropyl alcohol 70%	53315	1
Hedinger	Isopropyl alcohol 70%	53315	4
Hedinger	Isopropyl alcohol 70%	53316	1
AHD/Hedinger	Isopropyl alcohol 70%	53318	1
Hedinger	Isopropyl alcohol 70%	53320	1
Hedinger / Noweda	Isopropyl alcohol 70%	53321	1
Hedinger/Geilenkirchen	Isopropyl alcohol 70%	53320	1
Hedinger/Noweda	Isopropyl alcohol 70%	53320	1
Hedinger/AHD	Isopropyl alcohol 70%	53322	1
Hedinger/Gehe	Isopropyl alcohol 70%	53322	1
Hedinger/Alliance	Isopropyl alcohol 70%	53324	1
Hedinger / Noweda	Isopropyl alcohol 70%	53324	1
Hedinger	Isopropyl alcohol 70%	53324	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53324	1
Hedinger GmbH & Co. KG/...	Isopropyl alcohol 70%	53328	1
Hedinger/Fiebig	Isopropyl alcohol 70%	53326	1
Phoenix/AHD	Isopropyl alcohol 70%	53329	1
Hedinger / Phönix	Isopropyl alcohol 70%	53326	1
Hedinger GmbH & Co. KG/...	Isopropyl alcohol 70%	53329	1
AHD/Hedinger	Isopropyl alcohol 70%	53329	1
Caelo	Isopropyl alcohol 70%	53332	1
Hedinger/Phoenix	Isopropyl alcohol 70%	53332	1
Hedinger	Isopropyl alcohol 70%	53333	1
Hedinger/Gehe	Isopropyl alcohol 70%	53333	1
Noweda	Isopropyl alcohol 70%	53333	1
Hedinger/Geilenkirchen	Isopropyl alcohol 70%	53333	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Isopropyl alcohol 70%	53334	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53337	1
HEDINGER / Noweda	Isopropyl alcohol 70%	53340	1
Hedinger	Isopropyl alcohol 70%	53340	3
Hedinger / Gehe	Isopropyl alcohol 70%	53337	1
Hedinger/AHD	Isopropyl alcohol 70%	53340	1
Hedinger/Gehe	Isopropyl alcohol 70%	53338	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53341	1
Hedinger /Sanacorp WE: ...	Isopropyl alcohol 70%	53338	1
Hedinger/Gehe	Isopropyl alcohol 70%	53340	1
Phönix/Hedinger	Isopropyl alcohol 70%	53341	1
Fiebig/ Hedinger	Isopropyl alcohol 70%	53343	1
Hedinger/Gehe	Isopropyl alcohol 70%	53346	1
Hedinger/AHD	Isopropyl alcohol 70%	53346	1
Hedinger/Noweda	Isopropyl alcohol 70%	53343	1
Hedinger/Phönix	Isopropyl alcohol 70%	53348	1
Hedinger/Gehe	Isopropyl alcohol 70%	5335	1
Noweda/Hedinger	Isopropyl alcohol 70%	53350	1
Hedinger/Gehe	Isopropyl alcohol 70%	53352	1
Hedinger/Gehe	Isopropyl alcohol 70%	53351	1
Hedinger	Isopropyl alcohol 70%	53353	1
Hedinger	Isopropyl alcohol 70%	53354	3
Hedinger/Alliance	Isopropyl alcohol 70%	53353	1
Hedinger/ Sanacorp	Isopropyl alcohol 70%	53354	1
Hedinger/ Gehe	Isopropyl alcohol 70%	53354	1
Hedinger/AHD	Isopropyl alcohol 70%	53355	1
Hedinger/Noweda	Isopropyl alcohol 70%	53355	1
Hedinger/Gehe	Isopropyl alcohol 70%	53357	3
Hedinger/Noweda	Isopropyl alcohol 70%	53357	1
Hedinger / Kehr	Isopropyl alcohol 70%	53359	1
Hedinger	Isopropyl alcohol 70%	53359	2
Gehe/Hedinger26.04.16	Isopropyl alcohol 70%	53359	1
Hedinger/ Sanacorp	Isopropyl alcohol 70%	53360	1
Hedinger	Isopropyl alcohol 70%	53362	1
Hedinger	Isopropyl alcohol 70%	53363	1
Hedinger/Gehe	Isopropyl alcohol 70%	53367	1
Hedinger/Noweda	Isopropyl alcohol 70%	53360	1
Caelo	Isopropyl alcohol 70%	53367	1
Hedinger/Gehe	Isopropyl alcohol 70%	53366	1
Hedinger/Phoenix	Isopropyl alcohol 70%	53367	1
Hedinger/Noweda	Isopropyl alcohol 70%	53367	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53369	2
Hedinger	Isopropyl alcohol 70%	53368	2
Caelo	Isopropyl alcohol 70%	53369	1
Hedinger	Isopropyl alcohol 70%	53370	2
Hedinger	Isopropyl alcohol 70%	53371	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53370	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53371	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53372	1
Hedinger	Isopropyl alcohol 70%	53372	1
Gehe/Hedinger	Isopropyl alcohol 70%	53370	1
Gehe/Hedinger29.08.2016	Isopropyl alcohol 70%	53372	1
Hedinger	Isopropyl alcohol 70%	53373	1
Hedinger/ Kehr	Isopropyl alcohol 70%	53377	4
Hedinger/Sanacorp	Isopropyl alcohol 70%	53375	1
Hedinger/Gehe	Isopropyl alcohol 70%	53376	1
Lamotte/Sanacorp	Isopropyl alcohol 70%	53377	1
Hedinger/Phoenix	Isopropyl alcohol 70%	53377	1
Hedinger	Isopropyl alcohol 70%	53378	1
Gehe/Hedinger	Isopropyl alcohol 70%	53379	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger / AHD	Isopropyl alcohol 70%	53378	1
Caelo	Isopropyl alcohol 70%	53378	1
Hedinger/Gehe	Isopropyl alcohol 70%	53379	1
Hedinger/Noweda	Isopropyl alcohol 70%	53379	1
Hedinger/Noweda	Isopropyl alcohol 70%	53380	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53380	2
Hedinger/Gehe	Isopropyl alcohol 70%	53380	1
Geilenkirchen/Hedinger	Isopropyl alcohol 70%	53382	1
Hedinger/Phoenix	Isopropyl alcohol 70%	53382	1
Hedinger	Isopropyl alcohol 70%	53382	3
Hedinger	Isopropyl alcohol 70%	53383	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53382	1
Hedinger/Geibo	Isopropyl alcohol 70%	53383	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53383	1
Geilenkirchen/Hedinger	Isopropyl alcohol 70%	53383	1
Hedinger	Isopropyl alcohol 70%	53385	3
Hedinger/Noweda	Isopropyl alcohol 70%	53383	1
Caelo	Isopropyl alcohol 70%	53385	1
Hedinger / Noweda	Isopropyl alcohol 70%	53385	1
Hedinger/Phönix	Isopropyl alcohol 70%	53385	1
Aug. Hedinger	Isopropyl alcohol 70%	53385	2
Hedinger/Kehr	Isopropyl alcohol 70%	53386	1
Hedinger	Isopropyl alcohol 70%	53387	2
Hedinger/Gehe	Isopropyl alcohol 70%	53385	1
allianz	Isopropyl alcohol 70%	53388	1
WE: Phoenix 19.10.2016 ...	Isopropyl alcohol 70%	53388	2
HEDINGER/Phönix	Isopropyl alcohol 70%	53387	1
Hedinger	Isopropyl alcohol 70%	53388	1
Hedinger/Gehe	Isopropyl alcohol 70%	53388	1
Hedinger	Isopropyl alcohol 70%	53392	2
Hedinger	Isopropyl alcohol 70%	53390	1
Hedinger/Phönix	Isopropyl alcohol 70%	53392	1
Hedinger/Gehe	Isopropyl alcohol 70%	53392	1
Hedinger/Kehr	Isopropyl alcohol 70%	53393	1
Hedinger/Noweda	Isopropyl alcohol 70%	53395	1
Hedinger/Phönix	Isopropyl alcohol 70%	53397	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53396	1
Hedinger	Isopropyl alcohol 70%	53398	4
Sanacorp (Hedinger)	Isopropyl alcohol 70%	53397	1
Hedinger/Noweda	Isopropyl alcohol 70%	53397	1
hedinger/noweda	Isopropyl alcohol 70%	53397	1
Hedinger/ Sanacorp	Isopropyl alcohol 70%	53398	1
Caelo	Isopropyl alcohol 70%	53398	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53398	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53401	4
WE 27.12.16 Hedinger/Ph...	Isopropyl alcohol 70%	53401	1
Hedinger	Isopropyl alcohol 70%	53400	1
Hedinger/Gehe	Isopropyl alcohol 70%	53401	2
Hedinger/Phönix	Isopropyl alcohol 70%	53401	1
Hedinger	Isopropyl alcohol 70%	53401	1
Hedinger	Isopropyl alcohol 70%	53402	3
Hedinger	Isopropyl alcohol 70%	53404	3
Hedinger/Gehe	Isopropyl alcohol 70%	53402	2
Caelo	Isopropyl alcohol 70%	53401	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53404	1
Hedinger	Isopropyl alcohol 70%	53406	1
SanacorpWE:19.01.2017 ...	Isopropyl alcohol 70%	53407	2
Geilenkirchen/Hedinger	Isopropyl alcohol 70%	53404	1
Caelo	Isopropyl alcohol 70%	53408	1
Hedinger	Isopropyl alcohol 70%	53409	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
HEDINGER/Phönix	Isopropyl alcohol 70%	53410	1
Hedinger/Noweda	Isopropyl alcohol 70%	53410	1
Fagron	Isopropyl alcohol 70%	53411	3
Hedinger/Phoenix WE: 09. ...	Isopropyl alcohol 70%	53410	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53411	1
Hedinger	Isopropyl alcohol 70%	53411	1
Hedinger/ Phoenix	Isopropyl alcohol 70%	53413	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53412	1
Caelo	Isopropyl alcohol 70%	53416	1
Caelo	Isopropyl alcohol 70%	53413	1
Hedinger / Phönix	Isopropyl alcohol 70%	53413	1
Hedinger/Gehe	Isopropyl alcohol 70%	53416	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53416	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53416	3
Hedinger	Isopropyl alcohol 70%	53416	2
Hedinger	Isopropyl alcohol 70%	53421	1
hadinger/allianz	Isopropyl alcohol 70%	53420	1
Finteler Apotheke	Isopropyl alcohol 70%	53425	1
Hedinger	Isopropyl alcohol 70%	53419	1
Hedinger	Isopropyl alcohol 70%	53418	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53421	1
Caelo	Isopropyl alcohol 70%	53426	1
Hedinger/Gehe	Isopropyl alcohol 70%	53422	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53426	1
Franz GmbH, Phoenix	Isopropyl alcohol 70%	53426	1
Caelo	Isopropyl alcohol 70%	53428	1
Hedinger/Phönix	Isopropyl alcohol 70%	53426	1
Hedinger	Isopropyl alcohol 70%	53427	2
Hedinger/ Kehr	Isopropyl alcohol 70%	53429	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53432	4
Hedinger/Sanacorp	Isopropyl alcohol 70%	53432	1
Hedinger / Phönix	Isopropyl alcohol 70%	53430	1
Hedinger/Phönix	Isopropyl alcohol 70%	53432	1
Hedinger	Isopropyl alcohol 70%	53437	2
Caelo	Isopropyl alcohol 70%	53437	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53440	1
Hedinger	Isopropyl alcohol 70%	53439	1
Hedinger/Noweda	Isopropyl alcohol 70%	53439	1
Hedinger	Isopropyl alcohol 70%	53442	3
Hedinger/Kehr	Isopropyl alcohol 70%	53439	1
Hedinger	Isopropyl alcohol 70%	53443	1
Hedinger / Phönix	Isopropyl alcohol 70%	53442	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53443	1
Hedinger/ Noweda	Isopropyl alcohol 70%	53440	1
Hedinger/Noweda	Isopropyl alcohol 70%	53443	1
Hedinger	Isopropyl alcohol 70%	53444	3
Hedinger	Isopropyl alcohol 70%	53440	1
Hedinger/Phönix	Isopropyl alcohol 70%	53447	1
Hedinger	Isopropyl alcohol 70%	53449	2
Hedinger/Noweda	Isopropyl alcohol 70%	53454	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53451	2
Hedinger/GEHE	Isopropyl alcohol 70%	53453	1
Sanacorp WE: 25.10.17 1...	Isopropyl alcohol 70%	53452	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53456	1
HEDINGER/Phönix	Isopropyl alcohol 70%	53456	1
Hedinger/Phönix	Isopropyl alcohol 70%	53456	1
Hedinger	Isopropyl alcohol 70%	53458	2
Sanacorp WE: 21.11.17 1...	Isopropyl alcohol 70%	53460	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53460	1
Hedinger	Isopropyl alcohol 70%	53461	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger / Noweda	Isopropyl alcohol 70%	53461	1
Hedinger	Isopropyl alcohol 70%	53463	1
Hedinger	Isopropyl alcohol 70%	53465	1
Hedinger	Isopropyl alcohol 70%	53469	1
Hedinger	Isopropyl alcohol 70%	53462	1
Hedinger	Isopropyl alcohol 70%	53464	1
Hedinger	Isopropyl alcohol 70%	53470	2
Hedinger/Sanacorp	Isopropyl alcohol 70%	53470	1
Phönix E.K:5,68 -1L	Isopropyl alcohol 70%	53474	2
Hedinger/Gehe	Isopropyl alcohol 70%	53482	1
Hedinger	Isopropyl alcohol 70%	53480	1
Hedinger	Isopropyl alcohol 70%	53484	1
Hedinger	Isopropyl alcohol 70%	53491	1
Kehr/A.Hedinger	Isopropyl alcohol 70%	53490	1
Hedinger/Kehr	Isopropyl alcohol 70%	53497	1
Hedinger/Noweda	Isopropyl alcohol 70%	53485	1
Hedinger/Anzag	Isopropyl alcohol 70%	53497	1
Hedinger	Isopropyl alcohol 70%	53509	1
Caelo	Isopropyl alcohol 70%	6061401	1
Caelo	Isopropyl alcohol 70%	6081405	1
Eiegn/Defektur	Isopropyl alcohol 70%	616E-01080	1
Herstellung in Apotheke	Isopropyl alcohol 70%	70587B	1
Herstellung in Apotheke	Isopropyl alcohol 70%	70643B	1
Eigen/Defektur	Isopropyl alcohol 70%	716M-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	716M-01080	1
Gehe/Hedinger	Isopropyl alcohol 70%	83366	1
hochhaus-Apotheke	Isopropyl alcohol 70%	916I-01080	1
Sonnen - Apotheke	Isopropyl alcohol 70%	9/16	1
hochhaus-Apotheke	Isopropyl alcohol 70%	917M-01080	1
Labor	Isopropyl alcohol 70%		4
Rondell	Isopropyl alcohol 70%	1014M-07123	1
HirschApotheke	Isopropyl alcohol 70%	10151-01080	1
Caelo	Isopropyl alcohol 70%	1016A	1
Eigen Defektur	Isopropyl alcohol 70%	1016E-01080	1
Hubertus Apotheke	Isopropyl alcohol 70%	110416	1
Selbsterstellung	Isopropyl alcohol 70%	1017A-07096	1
Caelo	Isopropyl alcohol 70%	11031406	1
Hochhaus-Apotheke	Isopropyl alcohol 70%	1116E-01080	1
Caelo	Isopropyl alcohol 70%	11051501	1
Rondell Apotheke	Isopropyl alcohol 70%	1115Q-07123	1
Eigen/Defektur	Isopropyl alcohol 70%	1116M-01080	1
Rezeptur	Isopropyl alcohol 70%	1150617	1
Caelo	Isopropyl alcohol 70%	113363003	1
Hergestellt	Isopropyl alcohol 70%	1150727	1
Caelo	Isopropyl alcohol 70%	1160113	1
Defektur	Isopropyl alcohol 70%	12/15	1
Caelo	Isopropyl alcohol 70%	12031514	1
Caelo	Isopropyl alcohol 70%	12031515	1
Defektur	Isopropyl alcohol 70%	1214I	1
Defektur	Isopropyl alcohol 70%	1215A	1
Brinn- Holtz	Isopropyl alcohol 70%	1233-04/18	1
Hedinger	Isopropyl alcohol 70%		2
Wolf-Apotheke / Defektur	Isopropyl alcohol 70%	123456	2
Apotheke	Isopropyl alcohol 70%	12345	1
Caelo	Isopropyl alcohol 70%	123456789	1
sittig	Isopropyl alcohol 70%	130616	1
Rondell Apotheke	Isopropyl alcohol 70%	1316I-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	1316I-01080	2
Brinn-Holtz	Isopropyl alcohol 70%	1288-09/18	1
Caelo	Isopropyl alcohol 70%	13319404	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 70%	13363107	1
Gehe/Phönix	Isopropyl alcohol 70%	136363003	1
Defektur	Isopropyl alcohol 70%	14/14	1
Caelo	Isopropyl alcohol 70%	14026401	5
Caelo	Isopropyl alcohol 70%	14026402	5
Caelo	Isopropyl alcohol 70%	14026403	3
Caelo	Isopropyl alcohol 70%	140226402	1
Stern-Defektur	Isopropyl alcohol 70%	140310001	1
Hetterich	Isopropyl alcohol 70%	140403	1
Stern-Defektur	Isopropyl alcohol 70%	140404005	1
Caelo	Isopropyl alcohol 70%	14041501	1
STERN	Isopropyl alcohol 70%	140507002	1
Caelo	Isopropyl alcohol 70%	140505	1
STERN	Isopropyl alcohol 70%	140617002	1
Caelo	Isopropyl alcohol 70%	14076501	2
Caelo	Isopropyl alcohol 70%	14074205	1
Caelo	Isopropyl alcohol 70%	14076602	6
Caelo	Isopropyl alcohol 70%	14076603	2
Caelo	Isopropyl alcohol 70%	14101403	1
Caelo	Isopropyl alcohol 70%	14101402	1
Caelo	Isopropyl alcohol 70%	141117P801/17203002	1
Caelo	Isopropyl alcohol 70%	14114501	2
Hustadt Apotheke	Isopropyl alcohol 70%	141211	1
Caelo	Isopropyl alcohol 70%	14114601	2
Leyh Pharma	Isopropyl alcohol 70%	141311	1
Caelo	Isopropyl alcohol 70%	14146701	3
Caelo	Isopropyl alcohol 70%	14146703	1
Caelo	Isopropyl alcohol 70%	14146801	1
Caelo	Isopropyl alcohol 70%	14146801	3
Caelo	Isopropyl alcohol 70%	14146802	1
Caelo	Isopropyl alcohol 70%	14146803	6
Defektur	Isopropyl alcohol 70%	1414I	1
Defektur	Isopropyl alcohol 70%	1415M-01080	1
Eigenherstellung	Isopropyl alcohol 70%	1416E-01080	1
Caelo	Isopropyl alcohol 70%	1416702	1
Defektur/eigen	Isopropyl alcohol 70%	1416Q-01080	1
Caelo	Isopropyl alcohol 70%	141t46801	1
Caelo	Isopropyl alcohol 70%	14229101	4
Caelo	Isopropyl alcohol 70%	14229103	2
Sanacorp	Isopropyl alcohol 70%	14267805	1
Caelo	Isopropyl alcohol 70%	14330803	2
Caelo	Isopropyl alcohol 70%	14330802	1
Caelo	Isopropyl alcohol 70%	14330901	3
Caelo	Isopropyl alcohol 70%	14327203	1
Caelo	Isopropyl alcohol 70%	14330902	1
Caelo	Isopropyl alcohol 70%	14331003	3
Caelo	Isopropyl alcohol 70%	14331004	2
Leyh Pharma	Isopropyl alcohol 70%	143541	1
Caelo	Isopropyl alcohol 70%	15014102	3
Caelo	Isopropyl alcohol 70%	15014101	1
Caelo	Isopropyl alcohol 70%	15014204	1
Caelo	Isopropyl alcohol 70%	15014202	2
Caelo	Isopropyl alcohol 70%	15014301	1
Caelo	Isopropyl alcohol 70%	15014302	2
Sittig	Isopropyl alcohol 70%	150416	1
Caelo	Isopropyl alcohol 70%	15053505	1
Caelo	Isopropyl alcohol 70%	15076401	2
EuroPTC/Phönix	Isopropyl alcohol 70%	15076401	2
Fagron	Isopropyl alcohol 70%	15076501	1
Caelo	Isopropyl alcohol 70%	15076601	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 70%	15076602	7
Caelo	Isopropyl alcohol 70%	15076603	2
Schwanen-Apothek-Köster	Isopropyl alcohol 70%	151028-1	1
Hustadt-Apotheke	Isopropyl alcohol 70%	151204	1
Caelo	Isopropyl alcohol 70%	15150201	1
Defektur	Isopropyl alcohol 70%	1516M-01080	1
Caelo	Isopropyl alcohol 70%	15191902	2
Caelo	Isopropyl alcohol 70%	15192001	2
Caelo	Isopropyl alcohol 70%	15192002	2
Caelo	Isopropyl alcohol 70%	15192102	1
Caelo	Isopropyl alcohol 70%	15192101	2
Leyh Pharma	Isopropyl alcohol 70%	152121	1
Leyh Pharms	Isopropyl alcohol 70%	152521	1
Caelo	Isopropyl alcohol 70%	15248503	1
Caelo	Isopropyl alcohol 70%	15309001	4
Caelo	Isopropyl alcohol 70%	15309101	1
Caelo	Isopropyl alcohol 70%	15309201	2
Ley Pharma/Phönix	Isopropyl alcohol 70%	153741	1
Caelo	Isopropyl alcohol 70%	15403802	8
Caelo	Isopropyl alcohol 70%	15403902	1
Caelo	Isopropyl alcohol 70%	15430002	1
Caelo	Isopropyl alcohol 70%	15430004	1
Caelo	Isopropyl alcohol 70%	15437701	1
Caelo	Isopropyl alcohol 70%	15437801	1
Caelo	Isopropyl alcohol 70%	15437802	3
Leyh Pharma	Isopropyl alcohol 70%	154531	1
Leyh Pharma	Isopropyl alcohol 70%	155021	1
Schwanen-Apotheke, Köster	Isopropyl alcohol 70%	160104-1	1
Caelo	Isopropyl alcohol 70%	16049003	1
Caelo	Isopropyl alcohol 70%	16052014-100	1
Caelo	Isopropyl alcohol 70%	16071503	1
Caelo	Isopropyl alcohol 70%	16088201	2
Caelo	Isopropyl alcohol 70%	16088302	4
Caelo	Isopropyl alcohol 70%	16088403	2
Caelo	Isopropyl alcohol 70%	16121403	1
Caelo	Isopropyl alcohol 70%	16121402	1
Caelo	Isopropyl alcohol 70%	16138703	1
leyh pharma	Isopropyl alcohol 70%	161411	1
Hirsch Apotheke	Isopropyl alcohol 70%	1615I-01080	1
Caelo	Isopropyl alcohol 70%	16167601	2
Celo	Isopropyl alcohol 70%	16167802	1
Caelo	Isopropyl alcohol 70%	16167801	1
Ichthyol über Phoenix	Isopropyl alcohol 70%	16167802	1
Caelo	Isopropyl alcohol 70%	16167901	1
Defektur 1616E-01080	Isopropyl alcohol 70%	1616E-01080	1
Caelo	Isopropyl alcohol 70%	16202803	1
Caelo	Isopropyl alcohol 70%	16319001	1
Caelo	Isopropyl alcohol 70%	16319101	2
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	17/2015	1
Caelo	Isopropyl alcohol 70%	17082102	2
Caelo	Isopropyl alcohol 70%	17082302	1
Kastner	Isopropyl alcohol 70%	1714I-01080	1
Defektur	Isopropyl alcohol 70%	1716A-01080	1
Caelo	Isopropyl alcohol 70%	17358202	1
Caelo	Isopropyl alcohol 70%	17358201	1
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	18/16	1
Caelo	Isopropyl alcohol 70%	18031503	1
Caelo	Isopropyl alcohol 70%	18031504	1
Caelo	Isopropyl alcohol 70%	18061508	1
HirschApotheke	Isopropyl alcohol 70%	1815E-01080	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 70%	18061509	1
Caelo	Isopropyl alcohol 70%	18081510	1
Eigenherstellung	Isopropyl alcohol 70%	1815E-01080	1
Caelo	Isopropyl alcohol 70%	19111402	1
Caelo	Isopropyl alcohol 70%	18081511	1
Caelo	Isopropyl alcohol 70%	19111401	1
Eigenherstellung	Isopropyl alcohol 70%	1916A	1
Eigen/Defektur	Isopropyl alcohol 70%	1916M-01080	1
Eigenherstellung	Isopropyl alcohol 70%		3
Eigenherstellung	Isopropyl alcohol 70%	1953	1
Caelo	Isopropyl alcohol 70%	20011503	1
Caelo	Isopropyl alcohol 70%	2015181201E0	1
Caelo	Isopropyl alcohol 70%	20021506	1
Eigenherstellung	Isopropyl alcohol 70%	2015E-01080	1
Ottlik	Isopropyl alcohol 70%	2015090420	1
Rosen Apotheke	Isopropyl alcohol 70%	2017E-07125	1
Eigenherstellung	Isopropyl alcohol 70%	2016040520KS	1
Defektur/eigen	Isopropyl alcohol 70%	2216A-01080	1
Caelo	Isopropyl alcohol 70%	21114103	1
Eigenherstellung	Isopropyl alcohol 70%	2216E-01080	1
Defektur	Isopropyl alcohol 70%	2216B-01080	1
Defektur	Isopropyl alcohol 70%	22/14	1
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	23/15	2
Caelo	Isopropyl alcohol 70%	2216M-01080	1
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	23/16	1
Caelo	Isopropyl alcohol 70%	23061502	1
Caelo	Isopropyl alcohol 70%	23061501	1
Hubertus Apotheke	Isopropyl alcohol 70%	230915	1
Isopropylalkohol GMP/Gehe	Isopropyl alcohol 70%	2314M-01080	1
HirschApotheke	Isopropyl alcohol 70%	2315I-01080	1
Eigen/Defektur	Isopropyl alcohol 70%	2316I-01080	1
Frau Hennies	Isopropyl alcohol 70%	240216	1
Eigen	Isopropyl alcohol 70%	2416M-01080	1
hochhaus-apotheke/Eigen...	Isopropyl alcohol 70%	2417Q-01080	1
Caelo	Isopropyl alcohol 70%	25081401	1
Caelo	Isopropyl alcohol 70%	25091402	1
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600802T	1
Caelo	Isopropyl alcohol 70%	25091403	1
Eigen/Defektur	Isopropyl alcohol 70%	2516I-01080	1
noweda	Isopropyl alcohol 70%	2600787T	1
Winthrop AM/NOW	Isopropyl alcohol 70%	2600790T	1
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600830T	1
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600857T	1
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600807T	1
Winthrop	Isopropyl alcohol 70%	2600846T	1
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600846T	1
Winthrop/Sanacorp	Isopropyl alcohol 70%	2600882t	1
Winthrop/Sanacorp	Isopropyl alcohol 70%	2600873T	1
Winthrop/Phönix	Isopropyl alcohol 70%	2600852T	1
Zentiva/Alliance Health...	Isopropyl alcohol 70%	2600866T	1
Winthrop/Alliance Healt...	Isopropyl alcohol 70%	2600875T	1
Caelo	Isopropyl alcohol 70%	2600887T	1
Zentiva/Alliance Health...	Isopropyl alcohol 70%	2600887T	1
Winthrop/Alliance Healt...	Isopropyl alcohol 70%	2600892t	1
Sanacorp 26.11.15	Isopropyl alcohol 70%	2600905T	1
Zentiva/Alliance Health...	Isopropyl alcohol 70%	2600906T	2
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600915T	1
Winthrop Arzneimittel G...	Isopropyl alcohol 70%	2600916T	1
Caelo	Isopropyl alcohol 70%	2600907T	1
Sanacorp/Winthrop WE:18...	Isopropyl alcohol 70%	2600929T	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Winthrop	Isopropyl alcohol 70%	2600938T	1
Zentiva/Alliance Health...	Isopropyl alcohol 70%	2600948T	1
Winthrop	Isopropyl alcohol 70%	2600977T	1
Caelo	Isopropyl alcohol 70%	26071401	1
Caelo	Isopropyl alcohol 70%	2601004T	1
Rondell Apotheke	Isopropyl alcohol 70%	2615E-07123	1
Zentiva/Alliance Health...	Isopropyl alcohol 70%	2600953T	1
Caelo	Isopropyl alcohol 70%	2616E-01080	1
Eigen/Defektur	Isopropyl alcohol 70%	2616Q-01080	1
Eigen/Defektur	Isopropyl alcohol 70%	2716A-01080	1
Caelo	Isopropyl alcohol 70%	2716I-01080	1
Stifts-Apotheke	Isopropyl alcohol 70%	2716Q-01080	1
Defektur	Isopropyl alcohol 70%	2717I-07125	1
Defektur	Isopropyl alcohol 70%	28/14	1
Eigenherstellung	Isopropyl alcohol 70%	28072014	1
Rondell Apotheke	Isopropyl alcohol 70%	2814E-07123	1
Apotheke	Isopropyl alcohol 70%	2815E-01080	1
Defektur/eigen	Isopropyl alcohol 70%	2816M-01080	1
Linden-Apotheke Görlitz	Isopropyl alcohol 70%	2816E-01080	1
M.Kastner	Isopropyl alcohol 70%	2914M-01080	1
Defektur	Isopropyl alcohol 70%	2915E	1
Caelo	Isopropyl alcohol 70%	2x053426	1
Caelo	Isopropyl alcohol 70%	2916A-01080	1
Linden-Apotheke Görlitz	Isopropyl alcohol 70%	2915Q-01080	1
Hedinger/Phönix	Isopropyl alcohol 70%	3030046	1
Eigen/Defektur	Isopropyl alcohol 70%	3016M-01080	1
Defektur/eigen	Isopropyl alcohol 70%	3016E-01080	1
Hedinger/Phönix	Isopropyl alcohol 70%	303035	1
Rondell Apotheke	Isopropyl alcohol 70%	3016E-01080	1
Hedinger/Phönix	Isopropyl alcohol 70%	303037	3
Hedinger/Phönix	Isopropyl alcohol 70%	303039	3
Hedinger/Phönix	Isopropyl alcohol 70%	303040	1
Hedinger/Phönix	Isopropyl alcohol 70%	303045	1
Hedinger/Phönix	Isopropyl alcohol 70%	303042	1
Caelo	Isopropyl alcohol 70%	3115E-01080	1
Eigenherstellung	Isopropyl alcohol 70%	3116I-01080	1
eigenherstellung	Isopropyl alcohol 70%	3116Q-03043	1
Eigenherstellung	Isopropyl alcohol 70%	3215E-01080	1
N.N.	Isopropyl alcohol 70%	3216E	2
Defektur/eigen	Isopropyl alcohol 70%	3216I-01080	1
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	33/15	1
M.Kastner	Isopropyl alcohol 70%	3314I-01080	1
Apotheke	Isopropyl alcohol 70%	3366	1
Defektur	Isopropyl alcohol 70%	34/14	1
Eigenherstellung	Isopropyl alcohol 70%	3414E-01080	1
Eigen/Defektur	Isopropyl alcohol 70%	3515A-01080	1
wir leben Apotheke Tespe	Isopropyl alcohol 70%	3415M-01080	1
Frau Morgulis	Isopropyl alcohol 70%	3614E-01080	1
Eigen/Defektur	Isopropyl alcohol 70%	3516M-01080	1
eigen	Isopropyl alcohol 70%	3615I-01080	1
Eigen	Isopropyl alcohol 70%	3615A-01080	1
Eigen	Isopropyl alcohol 70%	3616M-01080	1
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	37/15	1
Hochhaus-Apotheke	Isopropyl alcohol 70%	3715A-01080	1
Defektur	Isopropyl alcohol 70%	37/14	1
Defektur	Isopropyl alcohol 70%	3715Q	1
Eigene	Isopropyl alcohol 70%	3815A-01080	1
Rondell-Apotheke	Isopropyl alcohol 70%	3814Q-07123	1
Rondell Apotheke	Isopropyl alcohol 70%	3816A-01080	1
Labor	Isopropyl alcohol 70%	3866	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	40509/50Wo	1
Wolf-Apotheke / Defektur	Isopropyl alcohol 70%	40729/50Hk	1
Kastner (eigen)	Isopropyl alcohol 70%	3914A-01080	1
Eigen	Isopropyl alcohol 70%	3914A-01080	1
Fischar	Isopropyl alcohol 70%	4114041	1
Wolf-Apotheke / Defektur	Isopropyl alcohol 70%	41125/1Hk	1
Otto Fischar GmbH & Co...	Isopropyl alcohol 70%	4115012	1
Defektur	Isopropyl alcohol 70%	4114I	1
Fischar	Isopropyl alcohol 70%	4115062	1
Eigenherstellung	Isopropyl alcohol 70%	4114Q-01080	1
Fischar	Isopropyl alcohol 70%	4115061	1
Eigenherstellung	Isopropyl alcohol 70%	4115A-01080	1
Defektur	Isopropyl alcohol 70%	4115Q	1
Wolf-Apotheke / Defektur	Isopropyl alcohol 70%	41212/50Si	1
Eigen	Isopropyl alcohol 70%	4215E-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	4215E-07123	1
Eigen/Defektur	Isopropyl alcohol 70%	4215M-01080	1
Defektur	Isopropyl alcohol 70%	43/14	1
Defektur	Isopropyl alcohol 70%	4315Q-01080	1
Hochhaus-Apotheke	Isopropyl alcohol 70%	4315E-01080	1
Sonnen-Apotheke Wefensl...	Isopropyl alcohol 70%	44/15	1
Defektur 4416A-01080	Isopropyl alcohol 70%	4416A-01080	1
Eigenherstellung	Isopropyl alcohol 70%	4515Q-01080	1
Eigenrezeptur	Isopropyl alcohol 70%	4616A-01080	1
Defektur	Isopropyl alcohol 70%	4415A-01080	1
HirschApotheke	Isopropyl alcohol 70%	4514I-01080	1
Hirsch Apotheke	Isopropyl alcohol 70%	4714I-01080	1
hochhaus-Apotheke	Isopropyl alcohol 70%	4715E-01080	1
Eigenrezeptur	Isopropyl alcohol 70%	4616E-01080	1
Defektur	Isopropyl alcohol 70%	4714M	1
Eigen/Defektur	Isopropyl alcohol 70%	4715E-01080	1
Eigenherstellung	Isopropyl alcohol 70%	4715M	1
Alte Apotheke	Isopropyl alcohol 70%	4716E-01080	1
Eigenherstellung	Isopropyl alcohol 70%	4716Q-01080	1
Defektur 4814N-01080	Isopropyl alcohol 70%	4814N-01080	2
Defektur	Isopropyl alcohol 70%	4917A-01080	1
Apotärze	Isopropyl alcohol 70%	4815I-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	4815M-01080	1
Eigen/Defektur	Isopropyl alcohol 70%	4815Q-01080	1
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	50103/50Wlo	1
Eigen/Defektur	Isopropyl alcohol 70%	5015A-01080	1
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	50227/50Wlo	1
Eigenherstellung	Isopropyl alcohol 70%	5015Q-01080	1
Wolf-Apotheke / Defektur	Isopropyl alcohol 70%	50204/50Bz	1
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	50330/50Wlo	1
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	50313/50Wlo	1
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	50715/50Wlo	1
C.Lange	Isopropyl alcohol 70%	50817/50La	1
Eigenherstellung	Isopropyl alcohol 70%	5114M-01080	1
CL. LANGE	Isopropyl alcohol 70%	50921/51LA	1
Defektur / Wolf-Apotheke	Isopropyl alcohol 70%	50828/50Wlo	1
Defektur	Isopropyl alcohol 70%	50907/50Si	1
Defektur	Isopropyl alcohol 70%	5214A	1
Wolf- Apotheke / Defektur	Isopropyl alcohol 70%	50914/51Wlo	1
hochhaus-Apotheke	Isopropyl alcohol 70%	5315E-01080	1
Eigen	Isopropyl alcohol 70%	5315I-01080	1
Hedinger/AHD	Isopropyl alcohol 70%	53290	1
Defektur	Isopropyl alcohol 70%	60209/50Hk	1
C.Lange	Isopropyl alcohol 70%	60307/51La	1
Claudia Lange	Isopropyl alcohol 70%	60620/50La	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Claudia Lange	Isopropyl alcohol 70%	60620/51La	1
Wolf-Apotheke	Isopropyl alcohol 70%	60623/50Dre	1
Claudia Lange	Isopropyl alcohol 70%	60704/50La	1
Defektur/Wolf-Apotheke	Isopropyl alcohol 70%	61215/50Dre	1
Fagron	Isopropyl alcohol 70%		1
Apotheke	Isopropyl alcohol 70%	70214/51Dre	1
Defektur	Isopropyl alcohol 70%		5
Hofmann und Sommer	Isopropyl alcohol 70%	940817	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-030715-1D	1
Defektur	Isopropyl alcohol 70%	80625/50A1	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-150216-1D	1
Defektur	Isopropyl alcohol 70%	71016/52Dre	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-300714-1D	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-241115-1D	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-110716-2D	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-200117-1D	1
Caelo	Isopropyl alcohol 70%	AL313	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-191214-1D	1
eigene herstellung	Isopropyl alcohol 70%	AL314	2
Eigene Herstellung	Isopropyl alcohol 70%	AL315	1
Eigene Herstellung	Isopropyl alcohol 70%	AL318	1
Eigene Herstellung	Isopropyl alcohol 70%	AL319	1
Eigene Herstellung	Isopropyl alcohol 70%	AL316	1
Eigene Herstellung	Isopropyl alcohol 70%	AL317	1
Eigene Herstellung	Isopropyl alcohol 70%	AL320	1
eigene Herstellung	Isopropyl alcohol 70%	AL321	1
Hubertus Apotheke	Isopropyl alcohol 70%	ChB060115	1
Caelo	Isopropyl alcohol 70%	Ch14146701	1
Eigenherstellung	Isopropyl alcohol 70%	Eigenherstellung	2
Rezeptur	Isopropyl alcohol 70%	Rezeptur	1
selbst	Isopropyl alcohol 70%		2
Neugereut Apotheke	Isopropyl alcohol 70%		2
Caelo	Isopropyl alcohol 70%		4
Eigen	Isopropyl alcohol 70%		1
Eigenherstellung	Isopropyl alcohol 70%	9051212	1
Defektur	Isopropyl alcohol 70%	2/13	1
Defektur	Isopropyl alcohol 70%	2/14	1
Caelo	Isopropyl alcohol 70%	2081202	1
Defektur_JP_07012013	Isopropyl alcohol 70%	213A-07015	1
Caelo	Isopropyl alcohol 70%	3091302	1
Neugereut Apotheke	Isopropyl alcohol 70%	214Q-01080	1
ASS	Isopropyl alcohol 70%	314A01080	1
Rondell	Isopropyl alcohol 70%	414Q-07123	2
Defektur_JP_28012013	Isopropyl alcohol 70%	513A-07015	1
Defektur	Isopropyl alcohol 70%	514A	1
Gehe	Isopropyl alcohol 70%	53202	1
Noweda	Isopropyl alcohol 70%	5072012	1
Gehe	Isopropyl alcohol 70%	53215	1
Hedinger	Isopropyl alcohol 70%	53228	1
Hedinger	Isopropyl alcohol 70%	53226	1
Caelo	Isopropyl alcohol 70%	53241	2
Hedinger	Isopropyl alcohol 70%	53241	3
hedinger	Isopropyl alcohol 70%	53236	2
Hedinger	Isopropyl alcohol 70%	53242	2
Hedinger	Isopropyl alcohol 70%	53240	1
Hedinger	Isopropyl alcohol 70%	53250	1
Hedinger	Isopropyl alcohol 70%	53245	1
Euro OTC	Isopropyl alcohol 70%	53245	1
Hedinger	Isopropyl alcohol 70%	53253	2
Phönix	Isopropyl alcohol 70%	53252	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Isopropyl alcohol 70%	53258	1
Hedinger	Isopropyl alcohol 70%	53255	1
Hedinger	Isopropyl alcohol 70%	53259	2
Hedinger	Isopropyl alcohol 70%	53262	1
Hedinger/Geilenkirchen	Isopropyl alcohol 70%	53262	1
Hedinger/AHD	Isopropyl alcohol 70%	53265	2
Hedinger	Isopropyl alcohol 70%	53267	2
Hedinger	Isopropyl alcohol 70%	53265	1
Hedinger	Isopropyl alcohol 70%	53268	2
Gehe	Isopropyl alcohol 70%	53268	1
Hedinger/Geilenkirchen	Isopropyl alcohol 70%	53266	1
Hedinger/AHD	Isopropyl alcohol 70%	53268	1
Hedinger	Isopropyl alcohol 70%	53274	1
Hedinger/Sanacorp	Isopropyl alcohol 70%	53300	1
Caelo	Isopropyl alcohol 70%	6041407	1
Hedinger/Hohe-Wart-Apot...	Isopropyl alcohol 70%	7082013/1	2
Defektur	Isopropyl alcohol 70%	80114Iso70	1
Caelo	Isopropyl alcohol 70%	7021307	1
Neugereut Apotheke	Isopropyl alcohol 70%	713M-01080	1
Defektur	Isopropyl alcohol 70%	7/14	1
Apoth.E.Bartel	Isopropyl alcohol 70%	70545B/2	1
Hirsch	Isopropyl alcohol 70%	714E-01080	1
Defektur	Isopropyl alcohol 70%	9/13	1
Doppeleiche Apotheke	Isopropyl alcohol 70%	913A-01080	1
Defektur	Isopropyl alcohol 70%	100513Iso70	1
Kastner	Isopropyl alcohol 70%	914I-01080	1
Defektur_JP_25022013	Isopropyl alcohol 70%	913A-07015	1
Hustadt Apotheke	Isopropyl alcohol 70%	121130	1
Apolda	Isopropyl alcohol 70%	120214	1
Caelo	Isopropyl alcohol 70%	121003	1
Herstellung	Isopropyl alcohol 70%	110416P2	1
Caelo	Isopropyl alcohol 70%	12121306	1
Hustadt Apotheke	Isopropyl alcohol 70%	120905	1
Klenk	Isopropyl alcohol 70%	1214M-1080	1
Caelo	Isopropyl alcohol 70%	12278601	7
Caelo	Isopropyl alcohol 70%	12101301	1
Caelo	Isopropyl alcohol 70%	12202801	2
Caelo	Isopropyl alcohol 70%	12136201	2
Caelo	Isopropyl alcohol 70%	12278704	4
Phönix	Isopropyl alcohol 70%	12338001	3
Caelo	Isopropyl alcohol 70%	12278702	2
Caelo	Isopropyl alcohol 70%	12338001	1
Caelo	Isopropyl alcohol 70%	12338104	4
Fagron	Isopropyl alcohol 70%	1234	2
Caelo	Isopropyl alcohol 70%	12338101	1
Caelo	Isopropyl alcohol 70%	13008902	2
Caelo	Isopropyl alcohol 70%	130009002	2
Fagron	Isopropyl alcohol 70%	13009002	2
Caelo	Isopropyl alcohol 70%	130009003	2
Caelo	Isopropyl alcohol 70%	13008903	2
Caelo	Isopropyl alcohol 70%	13009003	1
Caelo	Isopropyl alcohol 70%	130205	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	130206004	1
Stern-Defektur	Isopropyl alcohol 70%	130321004	1
Caelo	Isopropyl alcohol 70%	13031304	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	130327003	1
Stern-Defektur	Isopropyl alcohol 70%	130425004	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	130725001	1
Defektur_EW_04062013	Isopropyl alcohol 70%	130604002	1
Stern-Defektur	Isopropyl alcohol 70%	130805002	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Stern-Defektur	Isopropyl alcohol 70%	130703007	1
Defektur	Isopropyl alcohol 70%	130611/6	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	130813006	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	130912001	1
Caelo	Isopropyl alcohol 70%	13092201	1
Caelo	Isopropyl alcohol 70%	13092203	3
Phönix 12.06.2013	Isopropyl alcohol 70%	13092304	1
Stern-Defektur	Isopropyl alcohol 70%	130927002	1
Stern-Defektur	Isopropyl alcohol 70%	131108001	1
Stern-Rezeptur	Isopropyl alcohol 70%	131015001	1
Stern-Defektur	Isopropyl alcohol 70%	131216005	1
Caelo	Isopropyl alcohol 70%	13155401	2
Caelo	Isopropyl alcohol 70%	13155403	3
Unisan	Isopropyl alcohol 70%	13155502	1
STERN	Isopropyl alcohol 70%	131202002	1
Caelo	Isopropyl alcohol 70%	13155503	2
Ichthyol/Phönix	Isopropyl alcohol 70%	13155503	6
Caelo	Isopropyl alcohol 70%	13217602	1
Caelo	Isopropyl alcohol 70%	13217501	1
Caelo	Isopropyl alcohol 70%	13319405	2
Caelo	Isopropyl alcohol 70%	13217603	2
Caelo	Isopropyl alcohol 70%	13319502	2
Caelo	Isopropyl alcohol 70%	13337004	1
Caelo	Isopropyl alcohol 70%	13363003	2
Phönix 13.01.2014	Isopropyl alcohol 70%	13319404	1
Caelo	Isopropyl alcohol 70%	13363001	1
Stern-Defektur	Isopropyl alcohol 70%	140116004	1
Stern-Defektur	Isopropyl alcohol 70%	140210001	1
Caelo	Isopropyl alcohol 70%	14026303	2
Phonix	Isopropyl alcohol 70%	14026403	2
Caelo	Isopropyl alcohol 70%	14121205	1
Defektur	Isopropyl alcohol 70%	151113Iso70	1
Caelo	Isopropyl alcohol 70%	15111305	1
Caelo	Isopropyl alcohol 70%	15111304	1
Caelo	Isopropyl alcohol 70%	14101303	1
Defektur	Isopropyl alcohol 70%	16/13	1
Caelo	Isopropyl alcohol 70%	1613q-01080	1
Caelo	Isopropyl alcohol 70%	18061306	1
Defektur	Isopropyl alcohol 70%	180413Iso70	1
Purren Apotheke	Isopropyl alcohol 70%	18043	1
Caelo	Isopropyl alcohol 70%	18061307	1
hergestellt am 11.04.13	Isopropyl alcohol 70%	191012/925	1
Eigenherstellung	Isopropyl alcohol 70%	1823	2
Caelo	Isopropyl alcohol 70%	18101301	1
Hirsch	Isopropyl alcohol 70%	1913E-01080	1
Wensauer Apotheke	Isopropyl alcohol 70%		1
Caelo	Isopropyl alcohol 70%	20021414	1
Caelo	Isopropyl alcohol 70%	20021415	1
Defektur	Isopropyl alcohol 70%	200313Iso70	1
DEFEKTUR	Isopropyl alcohol 70%	2013/08/14-1	1
Caelo	Isopropyl alcohol 70%	20091302	1
Herstellung in Apotheke	Isopropyl alcohol 70%	20130613	1
Herstellung in Apotheke	Isopropyl alcohol 70%	20130102	1
Herstellung in Apotheke	Isopropyl alcohol 70%	20120911	1
Stern-Defektur	Isopropyl alcohol 70%	20130912	1
Stern-Defektur	Isopropyl alcohol 70%	20131022	1
DEFEKTUR	Isopropyl alcohol 70%	2013-09-16	1
Sonnen-Apotheke	Isopropyl alcohol 70%	210513H1	1
Herstellung in Apotheke	Isopropyl alcohol 70%	22042013	1
Caelo	Isopropyl alcohol 70%	22021303	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Defektur	Isopropyl alcohol 70%	210613Iso70	1
Caelo	Isopropyl alcohol 70%	22051302	1
Caelo	Isopropyl alcohol 70%	22051301	1
Caelo	Isopropyl alcohol 70%	221111206	1
Defektur	Isopropyl alcohol 70%	22/13	1
Caelo	Isopropyl alcohol 70%	221111205	1
geiger	Isopropyl alcohol 70%	2213E-01080	1
Caelo	Isopropyl alcohol 70%	23041304	1
Caelo	Isopropyl alcohol 70%	23041305	1
Defektur	Isopropyl alcohol 70%	23712	1
Caelo	Isopropyl alcohol 70%	24011405	1
Caelo	Isopropyl alcohol 70%	24081207	1
Caelo	Isopropyl alcohol 70%	2512A-01080	1
Caelo	Isopropyl alcohol 70%	24081303	1
Caelo	Isopropyl alcohol 70%	24101203	1
Winthrop	Isopropyl alcohol 70%	2600670T	1
Defektur	Isopropyl alcohol 70%	250713Iso70	1
Winthrop	Isopropyl alcohol 70%	2600685T	1
Winthrop	Isopropyl alcohol 70%	2600681T	1
Caelo	Isopropyl alcohol 70%	2600714T	1
GUMO/04.04.2013	Isopropyl alcohol 70%	2600726T	1
Eigenbedarf (Hergestel...	Isopropyl alcohol 70%	260613	1
Rondell Apotheke	Isopropyl alcohol 70%	2613Q-01080	2
Bombastus	Isopropyl alcohol 70%	2813Q-01080	1
Defektur	Isopropyl alcohol 70%	27/13	1
Defektur	Isopropyl alcohol 70%	280513Iso70	1
Defektur	Isopropyl alcohol 70%	2812M-01080	1
Defektur	Isopropyl alcohol 70%	290813Iso70	1
Doppeleiche Apotheke	Isopropyl alcohol 70%	3013A-01080	1
eigene Herstellung	Isopropyl alcohol 70%	3013M-01080	1
Hedinger/Anzag	Isopropyl alcohol 70%	303006	1
Hedinger	Isopropyl alcohol 70%	303007	2
Germania Apo	Isopropyl alcohol 70%	315032	1
Doppeleiche Apotheke	Isopropyl alcohol 70%	3413M-01080	1
Rondell Apotheke	Isopropyl alcohol 70%	3113M-07123	2
rondell-Apotheke	Isopropyl alcohol 70%	3413M-07123	2
Gutenberg-Apotheke	Isopropyl alcohol 70%	3212I-01080	1
Neugereut Apotheke	Isopropyl alcohol 70%	3531Q-01080	1
Defektur	Isopropyl alcohol 70%	37/13	1
Igel Apotheke	Isopropyl alcohol 70%	3813E-01080	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	3812E-01080	1
Defektur	Isopropyl alcohol 70%	3812	1
Apotheke am Holzweg	Isopropyl alcohol 70%	3912Y-01080	1
Doppeleiche Apotheke	Isopropyl alcohol 70%	4012M-01080	1
Euro OTC	Isopropyl alcohol 70%	4012Q-01080	1
Hirsch	Isopropyl alcohol 70%	4213E-01080	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	4212A-01080	1
Defektur	Isopropyl alcohol 70%	42/13	1
Hirsch	Isopropyl alcohol 70%	4413I-01080	1
Caelo	Isopropyl alcohol 70%	4813E-01080	1
Gutenberg-Apotheke	Isopropyl alcohol 70%	4612I-01080	1
Eigen	Isopropyl alcohol 70%	4413A-01080	1
Doppeleiche Apotheke	Isopropyl alcohol 70%	4912M-01080	1
Neugereut Apotheke	Isopropyl alcohol 70%	4913M-01080	1
Defektur	Isopropyl alcohol 70%	5013A	1
Hirsch	Isopropyl alcohol 70%	5013Q-01080	1
Defektur	Isopropyl alcohol 70%	53/12	1
Defektur	Isopropyl alcohol 70%	51/13	1
Sieg-Apotheke	Isopropyl alcohol 70%	661302027	1
Defektur	Isopropyl alcohol 70%	58/13	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl alcohol 70%	7581	1
Labor	Isopropyl alcohol 70%	7148	1
Apotheke	Isopropyl alcohol 70%	6710	1
Apotheke	Isopropyl alcohol 70%	8375	1
Caelo	Isopropyl alcohol 70%	8561	1
Labor	Isopropyl alcohol 70%	7988	1
Apotheke	Isopropyl alcohol 70%	9385	1
Apotheke	Isopropyl alcohol 70%	8441	1
Labor	Isopropyl alcohol 70%	8821	1
Apotheke	Isopropyl alcohol 70%	9993	1
Labor	Isopropyl alcohol 70%	8885	1
Martinus-Apotheke	Isopropyl alcohol 70%	A-171013-1D	2
Martinus-Apotheke	Isopropyl alcohol 70%	A-150314-1D	2
Martinus-Apotheke	Isopropyl alcohol 70%	A-131213-1D	2
eigene herstellung	Isopropyl alcohol 70%	a1309	1
eigene herstellung	Isopropyl alcohol 70%	a1308	1
Eigene Herstellung	Isopropyl alcohol 70%	AL311	1
Eigene Herstellung	Isopropyl alcohol 70%	AL312	1
Eigenherstellung	Isopropyl alcohol 70%	eigene	2
Hirsch	Isopropyl alcohol 70%	Die16041301	1
Caolo	Isopropyl alcohol 70%	D07104T1	1
Gehe/Wintrop	Isopropyl alcohol 70%	Ch2600715T	1
Hedinger	Isopropyl alcohol 70%	B230715	1
Löwen Apotheke	Isopropyl alcohol 70%	H28122012P1	1
Caelo	Isopropyl alcohol 70%	H170315T1	1
Caelo	Isopropyl alcohol 70%	Eigenherstellung	1
Caelo	Isopropyl alcohol 70%	H030315T1	1
Defektur	Isopropyl alcohol 70%	H17042013B1	1
Caelo	Isopropyl alcohol 70%	H22062013H1	1
MAROS Arznei, Noweda	Isopropyl alcohol 70%		1
Hedinger / Noweda	Isopropyl alcohol 70%	53266	1
Hedinger/Noweda	Isopropyl alcohol 70%	53266	1
Apotheke	Isopropyl alcohol 70%	6062013	1
Hedinger / Noweda	Isopropyl alcohol 70%	53267	1
Hubertus Apotheke	Isopropyl alcohol 70%	71014	1
Wintrop	Isopropyl alcohol 70%	80616-12	1
Sonnenberg-Apotheke	Isopropyl alcohol 70%	100214	1
Phönix	Isopropyl alcohol 70%	1217603	1
Hedinger	Isopropyl alcohol 70%	1313E-01079	1
Leyh Pharma	Isopropyl alcohol 70%	132021	1
Phönix	Isopropyl alcohol 70%	13217601	1
Caelo	Isopropyl alcohol 70%	13363104	1
Hetterich	Isopropyl alcohol 70%	140301	1
Defektur	Isopropyl alcohol 70%	1467/2014072920	1
Euro OCT	Isopropyl alcohol 70%	200214	1
Winthrop	Isopropyl alcohol 70%	2600706T	1
Wintrop AM GmbH/Anzag	Isopropyl alcohol 70%	2600754T	1
winthrop	Isopropyl alcohol 70%	2600735T	1
Wintrop	Isopropyl alcohol 70%	2600908t	1
Wintrop	Isopropyl alcohol 70%	2600869t	1
Wintrop	Isopropyl alcohol 70%	2600863T	2
Wintrop	Isopropyl alcohol 70%	2600969T	1
Defektur	Isopropyl alcohol 70%	3713Q-01080	1
Fischar	Isopropyl alcohol 70%	4116044	1
Otto Fischar gmbh	Isopropyl alcohol 70%	4116045	1
Apotheke	Isopropyl alcohol 70%	4813Q-0180	1

- 89 973 spectra from 501 *Apo-Ident* customers from a total of 26 162 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Isopropyl alcohol 70%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Isopropyl alcohol 70%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	760	0	93 114
Type C	0	1330	49	89 973

The substance/substance group *Isopropyl alcohol 70%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 99.2105 %)
Type C	100.0000 % (> 99.6586 %)	96.4467 % (> 96.2292 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31388	31388	0.00	76.05
31622	33105	6.00	67.01
31831	31388	3.94	73.38
33105	33105	0.00	66.96

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Isopropyl myristate
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30977-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Isopropyl myristate; Isopropylis myristas; Propylum myristinicum

Special notes

When selecting the *Isopropyl myristate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Isopropyl myristate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Isopropyl myrist...	12171522	30977	80	1408077
Caelo	Isopropyl myrist...	13144311	31363	40	1412539
Euro OTC	Isopropyl myrist...	1303034-03	31530	40	20130502
Euro OTC	Isopropyl myrist...	1508035	31774	40	20150907*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 4 reference samples from the substance/substance group *Isopropyl myristate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 344 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 207 spectra from 7 reference samples of the substance/substance group *Isopropyl myristate*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Isopropyl myristate	10119404	30193	7
Euro OTC	Isopropyl myristate	1303034-03	31530 [†]	20
Euro OTC	Isopropyl myristate	1508035	31774 [†]	20
Euro OTC	Isopropyl myristate	1608039	32776	40
Caelo	Isopropyl myristate	172729	33882	40
Euro OTC	Isopropyl myristate	1801035	33926	40
Caelo	Isopropyl myristate	181308	34234	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 667 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 200 spectra from 92 *Apo-Ident* customers from 57 batches from the substance/substance group *Isopropyl myristate*.
- These include spectra of independent samples from 53 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Euro OTC	Isopropyl myristate	1407015-03	13
Caelo	Isopropyl myristate	13144308	7
Euro OTC	Isopropyl myristate	1508035-01	13
Anzag	Isopropyl myristate	120100801	1
Euro OTC	Isopropyl myristate	1303034	1
Euro OTC	Isopropyl myristate	1303034-02	21
Euro/ Sanacorp	Isopropyl myristate	1303034-02	1
Euro OTG	Isopropyl myristate	1303034-02	1
Euro OTC	Isopropyl myristate	1303034-03	11
Euroo/Phoenix	Isopropyl myristate	1303034-03	1
Euro OTC	Isopropyl myristate	130303403	1
GEHE 26.11.2014	Isopropyl myristate	1303034-03	1
Fagron	Isopropyl myristate	13A17-B03-295774	2
Fagron	Isopropyl myristate	13A17-B05	1
Caelo	Isopropyl myristate	1407015-01	1
Euro OTC	Isopropyl myristate	1407015-01	4
Euro OTC	Isopropyl myristate	140701503	1
Noweda	Isopropyl myristate	1407015-01	1
Euro/ Sanacorp	Isopropyl myristate	1407015-01	1
Euro OTC	Isopropyl myristate	1407015-02	3
EuRho/Phönix	Isopropyl myristate	1407015-03	1
Phoenix	Isopropyl myristate	1407015-03	2
Caelo	Isopropyl myristate	14155902	1
Caelo	Isopropyl myristate	14155909	7
Caelo	Isopropyl myristate	14155914	2
EuRho/Noweda	Isopropyl myristate	1417015-03	1
Caelo	Isopropyl myristate	15035906	1
Euro/Phönix	Isopropyl myristate	1508035	1
Euro OTC	Isopropyl myristate	1508035	1
Caelo	Isopropyl myristate	1508035-01	2
Wolf-Apotheke	Isopropyl myristate	1508035-02	1
Caelo	Isopropyl myristate	15248304	2
	Isopropyl myristate	1508035-02	1
Euro OTC	Isopropyl myristate	1508035-02	1
Caelo	Isopropyl myristate	15248306	2
Caelo	Isopropyl myristate	15284804	1
Caelo	Isopropyl myristate	15284807	1
Caelo	Isopropyl myristate	15438204	3
Caelo	Isopropyl myristate	15438024	1
Euro OTC	Isopropyl myristate	1608039	1
Caelo	Isopropyl myristate	15438215	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Isopropyl myristate	1608039-01	10
Caelo	Isopropyl myristate	15438206	1
Euro OTC	Isopropyl myristate	1608039-02	1
Caelo	Isopropyl myristate	17272901	1
Caelo	Isopropyl myristate	30051502	1
Caelo	Isopropyl myristate	17359501	1
Euro OTC	Isopropyl myristate	3515Q-02898	1
Euro OTC	Isopropyl myristate	41140874	2
Euro OTC	Isopropyl myristate	1106026-01	1
Noweda	Isopropyl myristate	1106026-01	2
Euro OTC	Isopropyl myristate	968099	2
Caelo	Isopropyl myristate	11178408	1
EuRho	Isopropyl myristate	1106026-01	1
Caelo	Isopropyl myristate	11178415	1
Noweda	Isopropyl myristate	120100801	1
Caesar&Loretz GmbH D-40...	Isopropyl myristate	1201008-01	1
Noweda	Isopropyl myristate	1201008-01	4
Euro OTC	Isopropyl myristate	1201008-01	3
Fiebig ek: 3,00EUR 31.08.12	Isopropyl myristate	1201008-01	1
Euro OTC	Isopropyl myristate	120100801	1
EuRho	Isopropyl myristate	1201008-01	1
EuRho	Isopropyl myristate	12100801	1
Noweda	Isopropyl myristate	1212031-01	4
Euro OTC	Isopropyl myristate	1212031-01	4
Anzag 1.03.13 ek:3,45EUR	Isopropyl myristate	121203101	2
EuRho/Phönix	Isopropyl myristate	1212031-01	1
Noweda	Isopropyl myristate	12171511	1
Sanacorp	Isopropyl myristate	12171511	1
Caelo	Isopropyl myristate	12171522	4
Caelo	Isopropyl myristate	12171518	1
Phönix	Isopropyl myristate	12171522	1
Caelo	Isopropyl myristate	12248600	1
Euro OTC	Isopropyl myristate	130133402	1
Caelo	Isopropyl myristate	12248601	1
Caelo	Isopropyl myristate	13019902	5
Caesar & Loretz GmbH	Isopropyl myristate	13019902	1
Caelo	Isopropyl myristate	13019910	3
Phönix	Isopropyl myristate	13019910	1
Euro OTC	Isopropyl myristate	1303034-01	5
AHD	Isopropyl myristate	1301992	1
Phönix, 18.07.2013	Isopropyl myristate	130303401	1
Caelo	Isopropyl myristate	1303034-01	1
Euro OTC	Isopropyl myristate	1407015	1

- 91 152 spectra from 501 *Apo-Ident* customers from a total of 26 953 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Isopropyl myristate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Isopropyl myristate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The

following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	207	0	93 667
Type C	0	197	3	91 152

The substance/substance group *Isopropyl myristate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.1014 %)
Type C	100.0000 % (> 99.6586 %)	98.5000 % (> 97.0000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30977	30977	0.00	30.53
31363	31363	0.00	29.56
31530	31774	1.41	32.65
31774	31774	0.00	33.06

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Isopropyl palmitate
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30868-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Isopropyl palmitate; Isopropylis palmitas; Propylum palmitinicum

Special notes

When selecting the *Isopropyl palmitate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Isopropyl palmitate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Isopropyl palmit...	12229109	30868	40	20120807
Bombastus	Isopropyl palmit...	284060	31268	40	20131119
Caelo	Isopropyl palmit...	14034806	31364	40	AR-15-FG-000372-01
Bombastus	Isopropyl palmit...	306390	33040	40	20161219

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Isopropyl palmitate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Isopropyl palmitate*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Isopropyl palmitate	16091706	32570	40
Bombastus	Isopropyl palmitate	306390	32995	40
Caelo	Isopropyl palmitate	180672	34079	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 273 spectra from 146 *Apo-Ident* customers from 71 batches from the substance/substance group *Isopropyl palmitate*.
- These include spectra of independent samples from 67 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Isopropyl palmitate	16091701	7
Caelo	Isopropyl palmitate	17021402	4
Bombastus	Isopropyl palmitate	291043	8
VDL;4,53EUR;23.06.15	Isopropyl palmitate	29314	1
Sanacorp, 4,85EUR, 14.1.17	Isopropyl palmitate	302669	1
Bombastus	Isopropyl palmitate	299989	11
Caelo	Isopropyl palmitate	60092401	1
Bombastus	Isopropyl palmitate	266038	3
Bombastus	Isopropyl palmitate	281766	6
Caelo	Isopropyl palmitate	11207309	2
GEHE	Isopropyl palmitate	269359	1
Bombastus	Isopropyl palmitate	108611	1
Caelo	Isopropyl palmitate	2081401	1
Caelo	Isopropyl palmitate	11207301	5
Caelo	Isopropyl palmitate	14034801	1
AHD	Isopropyl palmitate	14034801	1
Caelo	Isopropyl palmitate	14034806	5
Caelo	Isopropyl palmitate	14034807	11
Sanacorp	Isopropyl palmitate	14034807	1
Caelo	Isopropyl palmitate	14213111	2
Noweda	Isopropyl palmitate	14213111	1
Caelo	Isopropyl palmitate	15128803	9
Sanacorp	Isopropyl palmitate	15128803	1
Caelo	Isopropyl palmitate	15169701	11
Caerol/Phoenix	Isopropyl palmitate	15169701	1
Caelo	Isopropyl palmitate	15169702	2
Caelo	Isopropyl palmitate	15169705	2
Caelo	Isopropyl palmitate	16071504	1
Caelo	Isopropyl palmitate	16021504	1
Caesar/Noweda	Isopropyl palmitate	16091701	1
Sanacorp	Isopropyl palmitate	16091701	2
Gehe	Isopropyl palmitate	16091701	1
Caelo	Isopropyl palmitate	16091706	8
Caelo	Isopropyl palmitate	16333401	1
Phönix	Isopropyl palmitate	16333401	1
Caelo	Isopropyl palmitate	1633401	1
Bombastus	Isopropyl palmitate	17001786	1
Bionorica/Kehr	Isopropyl palmitate	17000844	1
Bombastus	Isopropyl palmitate	170703IP/306390	1
Caelo	Isopropyl palmitate	281766	2
Caelo	Isopropyl palmitate	18067201	1
Bombastus	Isopropyl palmitate	284060	7
Bombastus	Isopropyl palmitate	285604	5
Audor Pharma	Isopropyl palmitate	285604	1
Brückenapotheke	Isopropyl palmitate	285604	1
Bombastus	Isopropyl palmitate	288083	5
AHD	Isopropyl palmitate	288083	1
Caelo	Isopropyl palmitate	288083	3
Bombas über Sanacorp	Isopropyl palmitate	290363	1
Sanacorp	Isopropyl palmitate	290363	1
Bombastus	Isopropyl palmitate	290363	3
Caelo	Isopropyl palmitate	291043	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Isopropyl palmitate	293140	13
Caelo	Isopropyl palmitate	293140	1
bombast/Px	Isopropyl palmitate	293140	1
Fagron	Isopropyl palmitate	293140	1
Bombastus	Isopropyl palmitate	297284	11
Caelo	Isopropyl palmitate	297284	1
phönix	Isopropyl palmitate	299989	1
Bombastus	Isopropyl palmitate	3015I-02899	1
Ichthyol-Gesellschaft	Isopropyl palmitate	299989	1
Bombastus	Isopropyl palmitate	302669	11
Bombastus	Isopropyl palmitate	307838	5
Phönix	Isopropyl palmitate	309754	1
Bombastus	Isopropyl palmitate	4613Q-02899	1
Bombastus	Isopropyl palmitate	262190	1
Bombastus	Isopropyl palmitate	269359	2
Hedinger	Isopropyl palmitate	271265	1
Bombastus	Isopropyl palmitate	272307	3
Bombastus	Isopropyl palmitate	271265	2
Caelo	Isopropyl palmitate	273823	2
Bombastus	Isopropyl palmitate	273823	2
Sanacorp	Isopropyl palmitate	273823	1
Caelo	Isopropyl palmitate	4061201	1
Caelo	Isopropyl palmitate	9324202	1
Bombastus	Isopropyl palmitate	9082013B	1
Fagron	Isopropyl palmitate	9G14-M08	2
Caelo	Isopropyl palmitate	12014106	2
Caelo	Isopropyl palmitate	12229109	2
Gehe	Isopropyl palmitate	12229112	1
Caelo	Isopropyl palmitate	131234905	1
Caelo	Isopropyl palmitate	13134903	4
Fiebig 21.08.13 EK5,60	Isopropyl palmitate	13134903	1
Caelo	Isopropyl palmitate	13134905	4
Noweda	Isopropyl palmitate	13134906	1
Caelo	Isopropyl palmitate	13134906	2
Caelo	Isopropyl palmitate	13134909	1
Euro OTC	Isopropyl palmitate		1
Caelo	Isopropyl palmitate	20210212-1	1
Caelo	Isopropyl palmitate	17091205	1
Caelo	Isopropyl palmitate	15071303	1
Caelo	Isopropyl palmitate	25031302	1
Bombastus	Isopropyl palmitate	26041301	1
Caelo	Isopropyl palmitate	26021402	1
Bombastus	Isopropyl palmitate	275014	7
Bombastus	Isopropyl palmitate	276178	2
Bombastus	Isopropyl palmitate	277740	3
Bombastus	Isopropyl palmitate	279468	4
Anzag	Isopropyl palmitate	277740Bombastus	1
Bombastus	Isopropyl palmitate	281057	6
Caelo	Isopropyl palmitate	28111309	1
Phoenix	Isopropyl palmitate	281766	1
Bombastus	Isopropyl palmitate	60092401	1
Gutenberg-Apotheke	Isopropyl palmitate	285604	1

- 91 079 spectra from 501 *Apo-Ident* customers from a total of 26 939 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Isopropyl palmitate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Isopropyl palmitate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	120	0	93 754
Type C	0	270	3	91 079

The substance/substance group *Isopropyl palmitate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6586 %)	98.9011 % (> 97.8022 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30868	31364	5.08	16.73
31268	31364	2.77	18.76
31364	31364	0.00	19.73
33040	31364	11.83	23.10

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Juniper oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30711-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Juniper oil; Juniperi aetheroleum; Oleum juniperi e baccararum

Special notes

When selecting the *Juniper oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Juniper oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Juniper oil	11379007	30711	40	1404599
Caelo	Juniper oil	13240811	31665	40	20130726
Caelo	Juniper oil	16160703	32918	40	1708747
Caelo	Juniper oil	171597	33311	40	20170601*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Juniper oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 70 spectra from 3 reference samples of the substance/substance group *Juniper oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Juniper oil	13240807	31244	10
Caelo	Juniper oil	13240811	31665 [†]	20
Caelo	Juniper oil	15185205	32352	40

- 93 804 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 20 spectra from 15 *Apo-Ident* customers from 12 batches from the substance/substance group *Juniper oil*.
- These include spectra of independent samples from 11 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
VDL;15.05.15;9,95EUR	Juniper oil	14115802	1
Caelo	Juniper oil	14115802	3
Caelo	Juniper oil	15185205	1
Caelo	Juniper oil	13240802	4
Caelo	Juniper oil	15141302	1
Caelo	Juniper oil	11379003	1
Caelo	Juniper oil	111275	1
Caelo	Juniper oil	11379006	1
Caelo	Juniper oil	11379012	2
Caelo	Juniper oil	11379007	1
Sanacorp	Juniper oil	11379012	1
Caelo	Juniper oil	12322707	1
Caelo	Juniper oil	11379017	1
Caelo	Juniper oil	18121306	1

- 91 332 spectra from 501 *Apo-Ident* customers from a total of 26 998 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Juniper oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Juniper oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	68	2	93 804
Type C	0	19	1	91 332

The substance/substance group *Juniper oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9762 %)	97.1429 % (> 92.8571 %)
Type C	100.0000 % (> 99.6590 %)	95.0000 % (> 80.0000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30711	30711	0.00	29.52
31665	30711	15.65	25.99
32918	32918	0.00	33.24
33311	33311	0.00	25.92

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Lactic acid**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30021-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Lactic acid; Acidum lacticum

Special notes

When selecting the *Lactic acid* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Lactic acid*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Lactic acid	12142108	30734	39	AR-14-FG-012230-01
Caelo	Lactic acid	13317311	31436	40	20131011
Caelo	Lactic acid	16144212	33271	40	20160527

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 119 spectra of 3 reference samples from the substance/substance group *Lactic acid*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 425 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Lactic acid*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Lactic acid	90646219	30021	40
Caelo	Lactic acid	13317311	31436 [†]	20
Caelo	Lactic acid	172176	33694	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 516 spectra from 273 *Apo-Ident* customers from 94 batches from the substance/substance group *Lactic acid*.
- These include spectra of independent samples from 91 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Lactic acid	11179209	1
Caelo	Lactic acid	11179217	7
Gutenberg-Apotheke	Lactic acid	13317303	1
Caelo	Lactic acid	13317308	9
Caelo	Lactic acid	13317311	17
Caelo	Lactic acid	13317313	5
Caelo	Lactic acid	14112112	5
Caelo	Lactic acid	13317315	5
Caelo	Lactic acid	14112115	10
A.H.D.	Lactic acid	15150402	1
Caelo	Lactic acid	15150410	15
Caelo	Lactic acid	15150401	6
AHD; 20.7.15; 5.84EUR	Lactic acid	15150402	1
Caelo	Lactic acid	15150413	6
Caelo	Lactic acid	15150421	11
Caelo	Lactic acid	15308608	18
Caelo	Lactic acid	12142102	8
Caelo	Lactic acid	16144201	17
Caelo	Lactic acid	13020201	19
Caelo	Lactic acid	12142112	12
Caelo	Lactic acid	13020203	7
Caelo	Lactic acid	13020211	19
Caelo	Lactic acid	13317301	32
GEHE	Lactic acid	21021180/2	1
Caelo	Lactic acid	11179201	2
Caelo	Lactic acid	11179207	4
Caelo	Lactic acid	12142162	1
Caelo	Lactic acid	13020212	1
Caelo	Lactic acid	13317307	13
Caelo	Lactic acid	13317303	2
Caesar & Loretz GmbH	Lactic acid	13317308	1
Sanacorp	Lactic acid	13317311	1
Cealo/Kehr	Lactic acid	13317312	1
Caelo	Lactic acid	14112017	1
Caelo	Lactic acid	13317393	1
Caelo	Lactic acid	14112106	3
Caelo	Lactic acid	14112107	12
Caesar & Loretz GmbH	Lactic acid	14112107	1
Noweda 18.05.2015	Lactic acid	141121119	1
Euro OTC	Lactic acid	14112112	1
Caesar & Loretz GmbH Hi...	Lactic acid	14112112	1
Fagron	Lactic acid	14112112	2
Caelo	Lactic acid	14112114	4
Phoenix Goettin	Lactic acid	14112115	1
Caelo	Lactic acid	14112116	3
Caelo	Lactic acid	14112119	17
Gehe	Lactic acid	14112119	1
Caelo	Lactic acid	1510402	1
Caelo	Lactic acid	1510421	1
Caelo	Lactic acid	15150402	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar&Loretz/Noweda	Lactic acid	15150402	1
Rondell Apotheke	Lactic acid	15150402	1
Caelo	Lactic acid	15150405	2
Caeolo/Noweda	Lactic acid	15150410	1
Noweda	Lactic acid	15150410	2
Caelo	Lactic acid	15150418	3
Caelo	Lactic acid	15308601	11
Caesar&Loretz GmbH,22.0...	Lactic acid	16144212	1
Caelo	Lactic acid	1062014C	1
Caelo	Lactic acid	17031502	1
Caelo	Lactic acid	2182	1
Caelo	Lactic acid	2101202	1
Caelo	Lactic acid	408Q-01034	1
Caelo	Lactic acid	5031301	1
Caelo	Lactic acid	9064607	1
Caelo	Lactic acid	913Q-01034	1
Caelo	Lactic acid	9321404	1
Caelo	Lactic acid	10/AKD1	1
Caelo	Lactic acid	10323202	1
Caelo	Lactic acid	11111644	1
Bombastus	Lactic acid	10323202	1
v.d.l	Lactic acid	11179207	1
Gehe	Lactic acid	11179207	1
Caelo	Lactic acid	11179214	1
Gehe	Lactic acid	11179217	1
Caelo	Lactic acid	12142101	10
Caelo	Lactic acid	12142108	7
Caelo	Lactic acid	12142104	2
Bombastus	Lactic acid	12142108	2
Caelo	Lactic acid	12142109	2
Phoenix	Lactic acid	1214212	1
AHD, 25.11.13	Lactic acid	1302011	1
Caesar/Lorenz	Lactic acid	13020201	1
Sanacorp	Lactic acid	13020201	1
Phönix, 21.03.2013	Lactic acid	13020201	1
Noweda	Lactic acid	13020208	1
Caelo	Lactic acid	13020208	7
Gehe	Lactic acid	13020208	1
Unisan	Lactic acid	13020208	1
Caelo	Lactic acid	1302021	1
Caesar&Loretz GmbH/Anzag	Lactic acid	13020211	1
Phönix	Lactic acid	13020211	1
Caelo	Lactic acid	13020808	1
Caelo	Lactic acid	130360	1
Caelo	Lactic acid	13030201	1
Caelo	Lactic acid	1309211034	1
Caelo	Lactic acid	13317304	3
Sanacorp	Lactic acid	13317307	1
Caelo	Lactic acid	13317309	2
Caesar&loretz,AHC	Lactic acid	13317309	1
Caelo	Lactic acid	133173082x	1
Euro OTC	Lactic acid	1503086-03	1
Caesar&Loretz GmbH	Lactic acid	15308601	1
Caelo	Lactic acid	1530860102263531	1
Caelo	Lactic acid	15308602	2
Phoenix	Lactic acid	15308608	1
Caelo	Lactic acid	15308609	1
Caelo	Lactic acid	15308613	11
Phoenix	Lactic acid	15308613	1
Rondell Apotheke	Lactic acid	15308613	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lactic acid	16121-01034	1
Caelo	Lactic acid	16141212	1
Sanacorp	Lactic acid	16144201	1
Gehe	Lactic acid	16144201	1
Caelo	Lactic acid	16144206	2
Caelo	Lactic acid	16144210	13
GEHE	Lactic acid	16144201	1
Caelo	Lactic acid	16144212	24
Noweda	Lactic acid	16144212	1
Caelo	Lactic acid	16144214	5
Caelo	Lactic acid	16144216	3
Caelo	Lactic acid	16144218	12
GEHE	Lactic acid	18034004	1
Caelo	Lactic acid	18034004	2
Caelo	Lactic acid	17217601	1
Caelo	Lactic acid	2014U-01034	1
Caelo	Lactic acid	20221110-4	1
Caelo	Lactic acid	2710Q-01034	1
	Lactic acid	71449138	2
Caelo	Lactic acid	81302029	1
Caelo	Lactic acid	90646369	1
Noweda	Lactic acid	81302139	1
Caelo	Lactic acid	A-220413-3	1
Caelo	Lactic acid	A-310513-1	1

- 90 836 spectra from 501 *Apo-Ident* customers from a total of 26 915 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Lactic acid* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Lactic acid* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	119	0	35 425
Type B	0	100	0	93 774
Type C	0	498	18	90 836

The substance/substance group *Lactic acid* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 94.9580 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6586 %)	96.5116 % (> 95.9302 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several

new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30734	30734	0.00	162.61
31436	30734	7.07	158.03
33271	30734	4.86	161.86

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Lanette® ointment (conserved)
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30499-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Lanette® ointment (conserved); Unguentum Lanette®

Special notes

When selecting the *Lanette® ointment (conserved)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Lanette® ointment (conserved)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Lanette® ointmen...	13257301	31166	40	1402587
Caelo	Lanette® ointmen...	13257301	31225	40	1403113
Caelo	Lanette® ointmen...	13393001	31454	30	20131206
Caelo	Lanette® ointmen...	153481	31913	40	20151111*
Caelo	Lanette® ointmen...	153482	31914	40	20151111*
Caelo	Lanette® ointmen...	160894	32093	40	20160318*
Caelo	Lanette® ointmen...	162580	32500	40	20160920*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 270 spectra of 7 reference samples from the substance/substance group *Lanette® ointment (conserved)*. These samples are listed above in the section *calibration samples*. The reference samples originate from 6 different batches.
- 35 274 spectra from a total of 834 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 870 spectra from 23 reference samples of the substance/substance group *Lanette® ointment (conserved)*.
- These include spectra of independent samples from 17 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Lanette® ointment (conserved)	11256001	30423	40
Caelo	Lanette® ointment (conserved)	90249059	30437	40
Caelo	Lanette® ointment (conserved)	11367801	30473	40
Caelo	Lanette® ointment (conserved)	11367901	30499	40
Caelo	Lanette® ointment (conserved)	12223801	30777	40
Caelo	Lanette® ointment (conserved)	13393001	31454†	10
Caelo	Lanette® ointment (conserved)	153481	31913†	20
Caelo	Lanette® ointment (conserved)	153482	31914†	20
Caelo	Lanette® ointment (conserved)	153483	31915	60
Caelo	Lanette® ointment (conserved)	162580	32501	40
Caelo	Lanette® ointment (conserved)	162581	32502	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Lanette [®] ointment (conserved)	162581	32503	40
Caelo	Lanette [®] ointment (conserved)	162582	32504	40
Caelo	Lanette [®] ointment (conserved)	162582	32505	40
Caelo	Lanette [®] ointment (conserved)	170567	33144	40
Caelo	Lanette [®] ointment (conserved)	170415	33145	40
Caelo	Lanette [®] ointment (conserved)	172430	33752	40
Caelo	Lanette [®] ointment (conserved)	172431	33753	40
Caelo	Lanette [®] ointment (conserved)	173110	34008	40
Caelo	Lanette [®] ointment (conserved)	180232	34196	40
Caelo	Lanette [®] ointment (conserved)	180870	34197	40
Caelo	Lanette [®] ointment (conserved)	181184	34198	40
Caelo	Lanette [®] ointment (conserved)	182334	34263	40

- 93 004 spectra from a total of 2146 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 453 spectra from 120 *Apo-Ident* customers from 103 batches from the substance/substance group *Lanette[®] ointment (conserved)*.
- These include spectra of independent samples from 100 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Lanette [®] ointment (conserved)	309918023301	1
GEHE	Lanette [®] ointment (conserved)	9639910	1
Caelo	Lanette [®] ointment (conserved)	13393001	17
Caelo	Lanette [®] ointment (conserved)	14200601	7
Caelo	Lanette [®] ointment (conserved)	14089801	14
Caelo	Lanette [®] ointment (conserved)	14277101	3
Caelo	Lanette [®] ointment (conserved)	14277301	5
Caelo	Lanette [®] ointment (conserved)	14277401	20
Caelo	Lanette [®] ointment (conserved)	15048801	1
Caelo	Lanette [®] ointment (conserved)	15048901	18
Caelo	Lanette [®] ointment (conserved)	15049001	12
Caelo	Lanette [®] ointment (conserved)	15095902	3
Caelo	Lanette [®] ointment (conserved)	15211201	6
Caelo	Lanette [®] ointment (conserved)	152113001	1
Caesar & Loretz GmbH / ...	Lanette [®] ointment (conserved)	15211301	1
Caelo	Lanette [®] ointment (conserved)	15211301	16
Caelo	Lanette [®] ointment (conserved)	15348201	2
Caelo	Lanette [®] ointment (conserved)	15211401	3
Caelo	Lanette [®] ointment (conserved)	152114021	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz GmbH	Lanette [®] ointment (conserved)	15348201	1
Caelo	Lanette [®] ointment (conserved)	15348301	8
Fagron	Lanette [®] ointment (conserved)	15348201	1
Sanacorp	Lanette [®] ointment (conserved)	160258001	1
Caelo	Lanette [®] ointment (conserved)	16089401	33
Caelo	Lanette [®] ointment (conserved)	160817UL	1
Sanacorp	Lanette [®] ointment (conserved)	16089401	1
Phoenix	Lanette [®] ointment (conserved)	16089402	2
Caelo	Lanette [®] ointment (conserved)	16089402	16
Caelo	Lanette [®] ointment (conserved)	16089404	2
Caelo	Lanette [®] ointment (conserved)	161104UL/16188701	1
GEHE	Lanette [®] ointment (conserved)	16184702	1
Caelo	Lanette [®] ointment (conserved)	16184702	6
Caelo	Lanette [®] ointment (conserved)	16188701	2
Caelo	Lanette [®] ointment (conserved)	16258001	15
Caelo	Lanette [®] ointment (conserved)	16188801	2
Gehe	Lanette [®] ointment (conserved)	16258001	1
Caelo	Lanette [®] ointment (conserved)	16258101	7
GEHE	Lanette [®] ointment (conserved)	16258101	2
Caelo	Lanette [®] ointment (conserved)	16258201	3
Caelo	Lanette [®] ointment (conserved)	17041501	5
Caelo	Lanette [®] ointment (conserved)	170328UL/16188801	1
Caelo	Lanette [®] ointment (conserved)	170530UL/16188801	1
Caelo	Lanette [®] ointment (conserved)	170710UL/16188801	1
Caelo	Lanette [®] ointment (conserved)	17134001	4
Caelo	Lanette [®] ointment (conserved)	17133901	2
Caelo	Lanette [®] ointment (conserved)	17133801	1
Caelo	Lanette [®] ointment (conserved)	171107UL/17134001	1
Caelo	Lanette [®] ointment (conserved)	17243001	1
Caelo	Lanette [®] ointment (conserved)	17243102	2
Caelo	Lanette [®] ointment (conserved)	17311001	2
Caelo	Lanette [®] ointment (conserved)	180102UL/17134001	1
Caelo	Lanette [®] ointment (conserved)	19082015C	1
Caelo	Lanette [®] ointment (conserved)	18023202	1
Caelo	Lanette [®] ointment (conserved)		1
Caelo	Lanette [®] ointment (conserved)	ChB280916	1
Caesar/Lorenz	Lanette [®] ointment (conserved)	12313102	1
Caelo	Lanette [®] ointment (conserved)	13046701	16
Caelo	Lanette [®] ointment (conserved)	13257301	14
Caelo	Lanette [®] ointment (conserved)	8041502	1
Caesar/Lorenz	Lanette [®] ointment (conserved)	13257501	1
Caelo	Lanette [®] ointment (conserved)	1339201	1
Caelo	Lanette [®] ointment (conserved)	13392701	11
Caelo	Lanette [®] ointment (conserved)	14089501	5
Phoenix	Lanette [®] ointment (conserved)	14089501	1
Caelo	Lanette [®] ointment (conserved)	14089601	4
AHD	Lanette [®] ointment (conserved)	14089801	1
Sanacorp	Lanette [®] ointment (conserved)	14089801	1
Caelo	Lanette [®] ointment (conserved)	14277201	3
Caelo	Lanette [®] ointment (conserved)	14221601	2
Caelo	Lanette [®] ointment (conserved)	1420061	1
Caelo	Lanette [®] ointment (conserved)	16071402	1
Caelo	Lanette [®] ointment (conserved)	20092014B	1
Caelo	Lanette [®] ointment (conserved)	280714B	1
Caelo	Lanette [®] ointment (conserved)	40/AKD2	1
Caelo	Lanette [®] ointment (conserved)	24062014A	1
Caelo	Lanette [®] ointment (conserved)	23012015C	1
Caelo	Lanette [®] ointment (conserved)	12081301	7
Caelo	Lanette [®] ointment (conserved)	11367901	4
Caelo	Lanette [®] ointment (conserved)	1208131	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lanette [®] ointment (conserved)	12223701	4
Caelo	Lanette [®] ointment (conserved)	1222380	1
Caelo	Lanette [®] ointment (conserved)	12223801	7
Caelo	Lanette [®] ointment (conserved)	1231302	1
Caelo	Lanette [®] ointment (conserved)	12313101	7
Caelo	Lanette [®] ointment (conserved)	13046601	5
GEHE	Lanette [®] ointment (conserved)	13046701	1
Caelo	Lanette [®] ointment (conserved)	13046802	5
Alliance	Lanette [®] ointment (conserved)	13046701	1
Caelo	Lanette [®] ointment (conserved)	13177501	12
Caelo	Lanette [®] ointment (conserved)	13234007	1
Caelo	Lanette [®] ointment (conserved)	13257501	13
Caelo	Lanette [®] ointment (conserved)	20021309	1
Caelo	Lanette [®] ointment (conserved)	20310812-5	1
Caelo	Lanette [®] ointment (conserved)	20110712-3	1
Caelo	Lanette [®] ointment (conserved)	28111311	1
Caelo	Lanette [®] ointment (conserved)	6012014D	1
Caelo	Lanette [®] ointment (conserved)	11367801	1
Caelo	Lanette [®] ointment (conserved)	613Q-03211	1
phönix, 20.11.2012	Lanette [®] ointment (conserved)	12081301	2
Phönix	Lanette [®] ointment (conserved)	12223701	1
Caelo	Lanette [®] ointment (conserved)	12170301	1
Caelo	Lanette [®] ointment (conserved)	12313001	3
Caelo	Lanette [®] ointment (conserved)	130456601	1
Caelo	Lanette [®] ointment (conserved)	1277	1
Caelo	Lanette [®] ointment (conserved)	12313102	2
Alliance Healthcare/Lam...	Lanette [®] ointment (conserved)	13046601	1
Caesar/Lorenz	Lanette [®] ointment (conserved)	13046701	1
Caelo	Lanette [®] ointment (conserved)	13177301	3
gehe	Lanette [®] ointment (conserved)	13177501	1
Caelo	Lanette [®] ointment (conserved)	13191201	6
Caesar/Lorenz	Lanette [®] ointment (conserved)	13257301	2
Caelo	Lanette [®] ointment (conserved)	13257402	1
Caelo	Lanette [®] ointment (conserved)	13392901	3
Caelo	Lanette [®] ointment (conserved)	1391201	1
Caelo	Lanette [®] ointment (conserved)	16102013C	1
Caelo	Lanette [®] ointment (conserved)	24022014C	1
Caelo	Lanette [®] ointment (conserved)	U-150513-3	1
Caelo	Lanette [®] ointment (conserved)	35151-03211	1
Anzah / Henry lammotte	Lanette [®] ointment (conserved)	13257501	1
Caelo	Lanette [®] ointment (conserved)	15211102	1
Caelo	Lanette [®] ointment (conserved)	15258001	1
Caelo	Lanette [®] ointment (conserved)	15348202	1
Caelo	Lanette [®] ointment (conserved)	12312101	1
Caelo	Lanette [®] ointment (conserved)	13046501	2
Caelo	Lanette [®] ointment (conserved)	12312901	1
Fagron	Lanette [®] ointment (conserved)	13046601	1

- 90 899 spectra from 501 *Apo-Ident* customers from a total of 26 910 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Lanette[®] ointment (conserved)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Lanette[®] ointment (conserved)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	270	0	35 274
Type B	0	869	1	93 004
Type C	0	397	56	90 899

The substance/substance group *Lanette[®] ointment (conserved)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.7778 %)
Type B	100.0000 % (> 99.9761 %)	99.8851 % (> 99.5402 %)
Type C	100.0000 % (> 99.6586 %)	87.6380 % (> 86.9757 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31166	31166	0.00	26.10
31225	31225	0.00	25.20
31454	31225	3.15	23.22
31913	31913	0.00	25.57
31914	31914	0.00	27.31
32093	32093	0.00	26.08
32500	32500	0.00	24.82

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Lanolin DAB**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30013-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Lanolin DAB; Lanolinum DAB; Unguentum adipis lanæ DAB

Special notes

When selecting the *Lanolin DAB* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Lanolin DAB*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Lanolin DAB	13062501	31003	40	20130226
Caelo	Lanolin DAB	14032001	31439	40	AR-15-FG-008283-01
Caelo	Lanolin DAB	15255201	31746	40	20150727
Fagron	Lanolin DAB	15G02-T11-024155	31794	40	AR-16-FG-004053-01
Bombastus	Lanolin DAB	294969	31795	40	20150611
Caelo	Lanolin DAB	15405502	31798	40	20151130
Caelo	Lanolin DAB	163026	32863	40	20161028*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 280 spectra of 7 reference samples from the substance/substance group *Lanolin DAB*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 7 different batches.
- 35 264 spectra from a total of 833 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 735 spectra from 20 reference samples of the substance/substance group *Lanolin DAB*.
- These include spectra of independent samples from 16 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Lanolin DAB	90870149	30013	30
Henry Lamotte	Lanolin DAB	4922600	30163	65
Bombastus	Lanolin DAB	00258524	30228	40
Henry Lamotte	Lanolin DAB	5264900	30922	40
Caelo	Lanolin DAB	15255201	31746 [†]	20
Fagron	Lanolin DAB	15602-T11-024155	31794 [†]	20
Bombastus	Lanolin DAB	294969	31795 [†]	20
Caelo	Lanolin DAB	15405502	31798 [†]	20
Caelo	Lanolin DAB	163022	32859	40
Caelo	Lanolin DAB	163023	32860	40
Caelo	Lanolin DAB	163024	32861	40
Caelo	Lanolin DAB	163025	32862	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Lanolin DAB	163444	33181	40
Caelo	Lanolin DAB	170184	33182	40
Caelo	Lanolin DAB	170187	33183	40
Caelo	Lanolin DAB	171351	33370	40
Caelo	Lanolin DAB	171354	33371	40
Caelo	Lanolin DAB	172965	34015	40
Caelo	Lanolin DAB	180443	34151	40
Caelo	Lanolin DAB	182302	34270	40

- 93 139 spectra from a total of 2147 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1589 spectra from 358 *Apo-Ident* customers from 299 batches from the substance/substance group *Lanolin DAB*.
- These include spectra of independent samples from 292 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Lanolin DAB	14165601	39
Caeolo/ Noweda	Lanolin DAB	14333701	1
Caelo	Lanolin DAB	12193702	2
Caelo	Lanolin DAB	12194201	4
Caelo	Lanolin DAB	14165901	29
Bombastus	Lanolin DAB	264721	1
Bombastus	Lanolin DAB	264723	2
Bombastus	Lanolin DAB	266983	1
Bombastus	Lanolin DAB	9/14	1
Bombastus	Lanolin DAB	270439	1
Caelo	Lanolin DAB	1042014A	1
Bombastus	Lanolin DAB	270443	1
Caelo	Lanolin DAB	1096249	2
Caelo	Lanolin DAB	213Q-01765	1
Caeol	Lanolin DAB	7022014A	1
Caelo	Lanolin DAB	7071205	1
Caelo	Lanolin DAB	713I-01765	1
Noweda	Lanolin DAB	7102013A	1
Caelo	Lanolin DAB	8072014A	1
Caelo	Lanolin DAB	8071503	1
Caelo	Lanolin DAB	9011406	1
Caelo	Lanolin DAB	8082014B	1
Caelo	Lanolin DAB	9112013B	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Eigene Herstellung	Lanolin DAB	10052016	1
Caelo	Lanolin DAB	9041401	1
Caelo	Lanolin DAB	11011404	1
Caelo	Lanolin DAB	10358301	1
Phönix	Lanolin DAB	11113902	1
Caelo	Lanolin DAB	11281701	2
DAB 2008	Lanolin DAB	11281701	1
Caelo	Lanolin DAB	11281702	2
Euro OTC	Lanolin DAB	11281901	2
Caelo	Lanolin DAB	11282001	3
Gehe	Lanolin DAB	11282001	1
Caelo	Lanolin DAB	11315101	4
Anzag	Lanolin DAB	11315101	1
Caelo	Lanolin DAB	11315102	3
Caelo	Lanolin DAB	12037002	1
Fagron	Lanolin DAB	11k18-t01	1
Phönix Gotha	Lanolin DAB	12037101	1
Caelo	Lanolin DAB	12037101	7
Caesar & Loretz GmbH, 4...	Lanolin DAB	12037101	1
VDL	Lanolin DAB	12037301	1
Caelo	Lanolin DAB	12037301	2
Caelo	Lanolin DAB	12037401	5
Caelo	Lanolin DAB	12065101	5
Noweda	Lanolin DAB	12065101	1
Caelo	Lanolin DAB	12065102	3
VDL	Lanolin DAB	12065101/32	1
Caelo	Lanolin DAB	12193701	23
Gehe	Lanolin DAB	12193701	1
Sanacorp	Lanolin DAB	12193701	2
Caesar & Loretz GmbH	Lanolin DAB	12193701	1
Caesar & Loretz GmbH	Lanolin DAB	12193701	1
Caelo	Lanolin DAB	12193901	12
Anzag/	Lanolin DAB	12193901	1
Phönix Gotha	Lanolin DAB	12193901	1
Großhandel	Lanolin DAB	12193901	1
Caelo	Lanolin DAB	12260301	5
Sanacorp	Lanolin DAB	12260301	1
Caesar & Loretz GmbH, 4...	Lanolin DAB	12260301	1
Caelo	Lanolin DAB	12260302	3
Caelo	Lanolin DAB	12308601	8
Caelo	Lanolin DAB	12308801	3
Sanacorp	Lanolin DAB	12308801	1
Caelo	Lanolin DAB	12308802	21
Phönix Gotha	Lanolin DAB	12308802	1
Sanacorp	Lanolin DAB	12308802	2
Erika apotheke	Lanolin DAB	12308802	1
Caelo	Lanolin DAB	123456	1
Bombastus	Lanolin DAB	124967	1
Caelo	Lanolin DAB	1230882	1
Fagron	Lanolin DAB	12C06-T02	1
Fagron	Lanolin DAB	12C06-T03	1
Fagron	Lanolin DAB	12E22-T04	5
Phönix 08.01.13	Lanolin DAB	12E22-T04	1
Fagron	Lanolin DAB	12E22-T06	1
Fagron	Lanolin DAB	12E22-T07	1
Kehr	Lanolin DAB	12H01-T01	1
Fagron	Lanolin DAB	12H01-T01	3
Caelo	Lanolin DAB	12H01-T01	1
Fagron	Lanolin DAB	12H01-T03	1
Caelo	Lanolin DAB	12H01-T02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lanolin DAB	12H01-T04	2
Sanacorp	Lanolin DAB	13062401	2
Phoenix	Lanolin DAB	12H01-T03	1
Caelo	Lanolin DAB	13062401	13
Phoenix 22.03.13	Lanolin DAB	13062401	2
Caesar & Loretz GmbH, 4...	Lanolin DAB	13062401	1
Phönix Gotha	Lanolin DAB	13062401	2
Alliance Healthcare	Lanolin DAB	13062401	1
Phönix	Lanolin DAB	13062401	1
Großhandel	Lanolin DAB	13062401	1
Caelo	Lanolin DAB	13062501	13
Fagron	Lanolin DAB	13062501	1
Gehe	Lanolin DAB	13062501	2
Caelo	Lanolin DAB	13062602	4
Caelo	Lanolin DAB	13062701	8
Caelo	Lanolin DAB	13062801	2
Caelo	Lanolin DAB	1306501	1
Caelo	Lanolin DAB	13102802	24
Alliance	Lanolin DAB	13102802	1
Noweda	Lanolin DAB	13102802	1
Caelo	Lanolin DAB	1310802	1
Bombastus	Lanolin DAB	13173601	1
phönix	Lanolin DAB	13121001	1
Caelo	Lanolin DAB	13227601	24
Gehe	Lanolin DAB	13227601	1
Caesar & Loretz GmbH	Lanolin DAB	13227601	1
Caelo	Lanolin DAB	13227602	5
Caelo	Lanolin DAB	13227701	17
Noweda	Lanolin DAB	13227701	1
Phönix Gotha	Lanolin DAB	13227701	2
Caelo	Lanolin DAB	13227801	7
Caelo	Lanolin DAB	13326501	27
Gehe am 21.2.14	Lanolin DAB	13326501	2
Caelo	Lanolin DAB	13326502	6
Caelo	Lanolin DAB	13326601	8
Gehe	Lanolin DAB	13326601	1
Krieger	Lanolin DAB	133266012x	1
Caelo	Lanolin DAB	13326701	30
Großhandel	Lanolin DAB	13326701	1
Phönix	Lanolin DAB	13326701	1
Caelo	Lanolin DAB	13345902	20
Caelo	Lanolin DAB	133326701	1
Caelo	Lanolin DAB	13345901	1
Noweda	Lanolin DAB	13345902	1
Fagron	Lanolin DAB	13B04-T20-00108	1
Fagron	Lanolin DAB	13804-T22-0002	1
Fagron	Lanolin DAB	13804T21001085	1
Fagron	Lanolin DAB	13B04-T20-001084	2
Fagron	Lanolin DAB	13B04-T21-001085	3
fargon	Lanolin DAB	13B04-T20-001084	1
Gehe	Lanolin DAB	13B04-T20-001084	1
phönix	Lanolin DAB	13B04-T20-001084	1
Fagron	Lanolin DAB	13B04-T20-0010B4	1
Fagron	Lanolin DAB	13B04-T22-001078	2
Fagron	Lanolin DAB	13B12-T01-001087	2
Fagron	Lanolin DAB	13E0760-0002	1
Fagron	Lanolin DAB	13E07-T18-003392	3
Caelo	Lanolin DAB	13E07-T20-003399	2
Fagron	Lanolin DAB	13E07-T21-003404	1
Fagron	Lanolin DAB	13E07-T20-003401	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lanolin DAB	13E07-T22-003407	1
Fagron	Lanolin DAB	13H13-12	1
Fagron	Lanolin DAB	13E07-Trg-oo3395	1
Fagron	Lanolin DAB	13H13-T08006173	1
Fagron	Lanolin DAB	13H13-T08-006173	2
	Lanolin DAB	13H13-T10-006179	1
Fagron	Lanolin DAB	13H13T11-006183	1
Fagron	Lanolin DAB	13H13T11006184	1
Fagron	Lanolin DAB	13H13-T12-006187	3
Fagron	Lanolin DAB	13H13-To9-006176	1
Caelo	Lanolin DAB	14011304	1
Caelo	Lanolin DAB	13H13-To8-006173	1
Caelo	Lanolin DAB	14032001	26
Caesar/Lorenz	Lanolin DAB	14032001	1
Phoenix 10.04.2014	Lanolin DAB	14032001	2
Gehe	Lanolin DAB	14032001	1
Caesar & Loretz GmbH	Lanolin DAB	14032001	1
sanacorp	Lanolin DAB	14032001	1
Gehe 29.07.14	Lanolin DAB	14032001	1
Caelo	Lanolin DAB	14032101	22
Großhandel	Lanolin DAB	14032101	1
Phönix Gotha	Lanolin DAB	14032101	1
Bombastus	Lanolin DAB	14032101	1
Caelo	Lanolin DAB	14032201	4
Noweda	Lanolin DAB	14032301	1
Caelo	Lanolin DAB	14032301	11
Sanacorp	Lanolin DAB	14032301	1
Caelo	Lanolin DAB	14032302	2
Caelo	Lanolin DAB	14032401	6
Caelo	Lanolin DAB	14101306	1
Bombastus	Lanolin DAB	14165601	1
Phönix	Lanolin DAB	14165601	1
Gehe	Lanolin DAB	14165601	2
Caesar & Loretz GmbH	Lanolin DAB	14165601	1
Caelo	Lanolin DAB	14165701	3
Noweda	Lanolin DAB	14165701	1
Sanacorp	Lanolin DAB	14165901	1
Phönix Gotha	Lanolin DAB	14165901	2
sanacorp	Lanolin DAB	14165901	1
Bombastus	Lanolin DAB	14166001	1
Caelo	Lanolin DAB	14166001	11
Fagron	Lanolin DAB	14166001	1
Caelo	Lanolin DAB	1433001	1
Caelo	Lanolin DAB	14333602	27
Sanacorp 31.12.14 EK 3,76	Lanolin DAB	14333602	2
Caesar&lorenz,Phoenix	Lanolin DAB	14333602	1
Sanacorp	Lanolin DAB	14333602	2
Caelo	Lanolin DAB	14333901	8
Caelo	Lanolin DAB	14333701	4
Noweda	Lanolin DAB	14333901	1
Caelo	Lanolin DAB	14334001	29
Phönix Gotha	Lanolin DAB	14334001	1
Sanacorp	Lanolin DAB	14334001	1
Caelo	Lanolin DAB	1433602	4
ANZAG 01.01.143	Lanolin DAB	14a06-t01-009011	1
Fagron	Lanolin DAB	14A06-T01-009011	4
Fagron	Lanolin DAB	14A06T01009011	1
Caelo	Lanolin DAB	14A06-T01-009011	1
Fagron	Lanolin DAB	14A06T02	3
Fagron	Lanolin DAB	14A06T02009014	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lanolin DAB	14A06-T02-009014	2
Fagron	Lanolin DAB	14A06-T03-009017	3
GEHE 26.11.2014	Lanolin DAB	14A06-T04-009021	1
Fagron	Lanolin DAB	14A06-T04-009021	1
Alliance health 12.01.15	Lanolin DAB	14A06-T05-0002	1
Fagron	Lanolin DAB	14A06-T06-009028	1
Fagron	Lanolin DAB	14E06T07	1
Fagron	Lanolin DAB	14E06-T07-012186	5
Fagron	Lanolin DAB	14E06-T07	1
Gehe	Lanolin DAB	14E06-T07-012186	1
Fagron	Lanolin DAB	14E06-T09-012194	5
Fagron	Lanolin DAB	14E06-T10	1
Fagron	Lanolin DAB	14K11-T01-016976	3
Fagron	Lanolin DAB	14E06-T10-012197	1
Caelo	Lanolin DAB	14K11-T01-016976	2
Fagron	Lanolin DAB	14E06-T11-012200	1
Fagron	Lanolin DAB	14K11-T03-016983	2
GEHE	Lanolin DAB	14K11-T03-016984	1
Fagron	Lanolin DAB	14K11-T03-016984	2
Caelo	Lanolin DAB	15024003	1
Caelo	Lanolin DAB	15024102	25
Phönix	Lanolin DAB	15024102	1
Caesar & Loretz GmbH	Lanolin DAB	15024102	1
Caelo	Lanolin DAB	15024201	28
Phönix Gotha	Lanolin DAB	15024201	2
Caelo	Lanolin DAB	15024301	2
Caelo	Lanolin DAB	15045502	1
Caelo	Lanolin DAB	1515Q-01765	1
Caelo	Lanolin DAB	15193401	18
Phönix	Lanolin DAB	15193401	1
Phönix Gotha	Lanolin DAB	15193401	1
Caelo	Lanolin DAB	15193501	24
Caeolo/Noweda	Lanolin DAB	15193501	1
Caelo	Lanolin DAB	15255201	20
Caelo	Lanolin DAB	1525501	1
Gehe	Lanolin DAB	15255201	1
Fagron	Lanolin DAB	15255201	1
Caelo	Lanolin DAB	15255202	22
Cael Sanacorp 05.02.2016	Lanolin DAB	15255202	1
Caelo	Lanolin DAB	15255301	6
Sanacorp	Lanolin DAB	15255301	1
anzag	Lanolin DAB	15255301	1
Castor/Phönix	Lanolin DAB	15255301	1
Noweda	Lanolin DAB	15255501	2
Caelo	Lanolin DAB	15255501	2
Caelo	Lanolin DAB	15255601	6
Noweda	Lanolin DAB	15255601	2
Caelo	Lanolin DAB	15405501	15
Alliance H.	Lanolin DAB	15405501	1
Caelo/ AHD	Lanolin DAB	15405501	1
Noweda	Lanolin DAB	15405501	1
Caelo	Lanolin DAB	15405502	26
caesar&lorenz,Phönix	Lanolin DAB	15405502	1
WE:18.03.2016 Sanacorp	Lanolin DAB	15405502	1
Caelo	Lanolin DAB	15405601	27
Fagron	Lanolin DAB	15405601	1
Caelo	Lanolin DAB	15405801	3
Caelo	Lanolin DAB	15c1903020815	1
Caelo	Lanolin DAB	15C19-T01-020809	2
Fagron	Lanolin DAB	15C19-T01-020809	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Lanolin DAB	15C19-T02-020812	1
Fargron/Noweda	Lanolin DAB	15C19-T01-020809	1
Fargron	Lanolin DAB	15C19-T02-020812	8
Fargron	Lanolin DAB	15c19-t02-020812	1
Fargron	Lanolin DAB	15C19-T03-020815	6
Fargron	Lanolin DAB	15C19T03	1
Fargron	Lanolin DAB	15c19-t03-020815	1
Fargron	Lanolin DAB	15c19-t05-020821	1
Fargron	Lanolin DAB	15E04-T01-022046	2
Fargron	Lanolin DAB	15G02-T10-024152	2
Fargron	Lanolin DAB	15C19-T04-020818	1
Fargron	Lanolin DAB	15C19-T05-020821	1
Fargron	Lanolin DAB	15C-9-T02020812	1
Caelo	Lanolin DAB	1602301	1
Caelo	Lanolin DAB	16011501	1
Caelo	Lanolin DAB	16054602	4
Caelo	Lanolin DAB	16054701	7
Gehe	Lanolin DAB	16054602	1
Ichthyol-Gesellschaft	Lanolin DAB	16054701	1
Caelo	Lanolin DAB	16054801	28
Caelo	Lanolin DAB	16054701	1
Caelo/AHD	Lanolin DAB	16054801	2
WE:06.04.2016 Phönix	Lanolin DAB	16054901	3
Caelo	Lanolin DAB	16054901	31
Phönix	Lanolin DAB	16054901	1
Phoenix	Lanolin DAB	16054901	2
Euro OTC	Lanolin DAB	16054901	2
caesar & Loretz	Lanolin DAB	16054901	1
Caelo	Lanolin DAB	16064901	1
phönix	Lanolin DAB	16189401	2
Caelo	Lanolin DAB	16189401	20
Caesar & Loretz	Lanolin DAB	16189401	1
Phoenix	Lanolin DAB	16189401	1
Caesar & Loretz GmbH Hi...	Lanolin DAB	16189401	1
SanacorpWE:12.12.2016 ...	Lanolin DAB	16189401	1
AHD	Lanolin DAB	16189401	1
Noweda	Lanolin DAB	16189401	1
Caelo	Lanolin DAB	16189402	1
Anzag	Lanolin DAB	16189501	1
Caelo	Lanolin DAB	16189501	20
GEHE	Lanolin DAB	16189501	2
Caesar & Loretz GmbH	Lanolin DAB	16189501	1
Caelo	Lanolin DAB	16189501/160314L	1
Alte Apotheke	Lanolin DAB	16189602	1
Caelo	Lanolin DAB	16189701	1
Caelo	Lanolin DAB	16302202	2
Caelo	Lanolin DAB	16302301	25
Phoenix	Lanolin DAB	16302301	1
SanacorpWE:27.01.2017 ...	Lanolin DAB	16302301	1
Caelo	Lanolin DAB	16302401	27
Caelo	Lanolin DAB	16302302	1
AHDE/Cael030.12.16 EUR5,33	Lanolin DAB	16302401	1
Caeolo/Noweda	Lanolin DAB	16302401	1
Hubertus Apotheke	Lanolin DAB	16302401	1
Bombastus	Lanolin DAB	16302401	1
Fargron	Lanolin DAB	16302401	1
Caelo	Lanolin DAB	16302501	1
Caelo	Lanolin DAB	16344401	11
Caesar & Loretz GmbH	Lanolin DAB	16344401	1
Fargron	Lanolin DAB	16G01-T03-034847	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lanolin DAB	16G01-T04-034850	3
Fagron	Lanolin DAB	16G01-T05-034853	7
Fagron	Lanolin DAB	16G01-T06-034857	1
Fagron	Lanolin DAB	16G01-T07-034861	1
Fagron	Lanolin DAB	16G01-T08-034867	4
Bombastus	Lanolin DAB	17001446	1
Bombastus	Lanolin DAB	17001447	1
26.07.2018	Lanolin DAB	17000364	1
Caelo	Lanolin DAB	17018401	29
Gehe	Lanolin DAB	17000364	1
Noweda	Lanolin DAB	17018401	1
GEHE	Lanolin DAB	17018401	1
Caesar & Loretz/Noweda	Lanolin DAB	17018401	1
Caelo	Lanolin DAB		1
Caelo	Lanolin DAB	17018701	12
GEHE	Lanolin DAB	17018701	1
Caesar & Loretz GmbH	Lanolin DAB	17018701	1
Caelo	Lanolin DAB	17082015D	1
Sanacorp	Lanolin DAB	17078401	1
Caelo	Lanolin DAB	17135101	13
phönix	Lanolin DAB	17135101	1
Phoenix WE:04.11.17 100...	Lanolin DAB	17135101	1
Sanacorp	Lanolin DAB	17135101	1
Noweda	Lanolin DAB	17135101	1
GEHE	Lanolin DAB	17135101	1
Caelo	Lanolin DAB	17135201	12
Cealo	Lanolin DAB	17222401	2
Caelo	Lanolin DAB	17222401	5
Caelo	Lanolin DAB	17222601	10
Phönix	Lanolin DAB	17222601	1
Sanacorp W.E: 31.03.20...	Lanolin DAB	17222601	1
Noweda	Lanolin DAB	17222601	1
Caesar & Loretz/ Noweda	Lanolin DAB	17222601	1
AHD	Lanolin DAB	17222601	1
Caelo	Lanolin DAB	17296501	1
Caelo	Lanolin DAB	17296601	1
Caelo	Lanolin DAB	17296602	3
Caelo	Lanolin DAB	1735201	1
Caelo	Lanolin DAB	17A03-T05-040010	1
Fagron	Lanolin DAB	17A03-T02-040010	1
Fagron	Lanolin DAB	17A03-T02-040001	1
Caelo	Lanolin DAB	17296702	1
Fagron	Lanolin DAB	17A03-T05-040010	2
Fagron	Lanolin DAB	17H24-T01-047688	1
Caelo	Lanolin DAB	18092014B	1
Caelo	Lanolin DAB	17H24-T04-047697	1
Fagron	Lanolin DAB	17A03-T03-040004	1
Fagron	Lanolin DAB	17H24T04047697	1
Caelo	Lanolin DAB	20240912-1	1
Caelo	Lanolin DAB	23012015B	1
Caelo	Lanolin DAB	23326501	1
Euro OTC	Lanolin DAB	260416v1	1
Caelo	Lanolin DAB	26061401	1
Caelo	Lanolin DAB	26052014B	1
Bombastus	Lanolin DAB	274414	1
Bombastus	Lanolin DAB	274412	1
Bombstus	Lanolin DAB	273293	1
Bombastus	Lanolin DAB	274413	2
AHD , 05.09.13	Lanolin DAB	274414	1
Bombastus	Lanolin DAB	276824	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Lanolin DAB	277967	2
Bombastus	Lanolin DAB	277969	1
Bombastus	Lanolin DAB	277970	4
Bombastus	Lanolin DAB	280111	3
AHD, 08.11.13	Lanolin DAB	280111	1
Bombastus	Lanolin DAB	282979	2
Raben-Apotheke Mildenau	Lanolin DAB	282976	1
Bombastus	Lanolin DAB	285096	2
Bombastus	Lanolin DAB	286688	1
Bombastus	Lanolin DAB	285097	4
Bombastus	Lanolin DAB	286690	1
Bombastus	Lanolin DAB	288634	2
Bombastus	Lanolin DAB	286691	3
Bombasatus	Lanolin DAB	288635	1
Bombastus	Lanolin DAB	288636	2
Bombastus	Lanolin DAB	288637	4
Caelo	Lanolin DAB	290316L	1
Bombastus	Lanolin DAB	290885	4
Bombastus	Lanolin DAB	292904	1
Bombastus	Lanolin DAB	292906	3
Bombastus	Lanolin DAB	294967	2
Bombastus	Lanolin DAB	297583	4
Bombastus	Lanolin DAB	299082	5
Bombastus	Lanolin DAB	300144	1
Bombastus	Lanolin DAB	301170	2
Bombastus	Lanolin DAB	302808	1
Bombastus	Lanolin DAB	302807	1
Bombastus	Lanolin DAB	304540	1
Bombastus	Lanolin DAB	302809	4
Bombastus	Lanolin DAB	304541	5
Euro OTC	Lanolin DAB	306635	1
Bombastus	Lanolin DAB	306634	1
Bombastus	Lanolin DAB	305440	1
Bombastus	Lanolin DAB	306635	1
GEHE	Lanolin DAB	309000	1
Bombastus	Lanolin DAB	309000	1
Caelo	Lanolin DAB	31121301	1
phönix	Lanolin DAB		1
Sanacorp	Lanolin DAB	34/AKD1	1
Caelo	Lanolin DAB	5113A-01765	1
Caelo	Lanolin DAB	3613A-01765	1
phönix	Lanolin DAB	34/AKD2	1
anzag, 05.11.2012	Lanolin DAB	5152000	1
Caelo	Lanolin DAB	4413M01765	1
Lamotte/Noweda	Lanolin DAB	5152000	1
Caelo	Lanolin DAB	5211Q-01765	1
Henry Lamotte Oils	Lanolin DAB	5264900	1
Anzag 10.04.13	Lanolin DAB	5264900	1
Henry Lamotte Oils GmbH...	Lanolin DAB	5264900	1
Anzag/Henry Lamotte	Lanolin DAB	5264900	1
Anzag	Lanolin DAB	5264900	1
Lamotte/Noweda	Lanolin DAB	5264900	1
Phönix	Lanolin DAB	5316400	1
Henry	Lanolin DAB	5316400	1
Henry Lamotte Oils / AHD	Lanolin DAB	5375500	3
Lamotte/Noweda	Lanolin DAB	5375500	1
Alliance Healthcare	Lanolin DAB	5375500	1
Henry Lamotte Oils GmbH...	Lanolin DAB	5375500	2
Henry Lamotte	Lanolin DAB	5375500	2
lamotte	Lanolin DAB	5447700	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte/Anzag	Lanolin DAB	5447700	1
HenryLamotte/Alliance	Lanolin DAB	5517200	1
Lamotte/Alliance	Lanolin DAB	5517200	1
Henry Lamotte/ Noweda	Lanolin DAB	5517200	1
AHD	Lanolin DAB	5517200	1
Henry Lamotte/ Noweda	Lanolin DAB	5651400	2
Henry/AHD	Lanolin DAB	5651400	1
Sanacorp/ Henry Lamotte...	Lanolin DAB	5651400	1
Lamotte/Alliance	Lanolin DAB	5651400	5
Lamotte	Lanolin DAB	5651400	1
Caelo	Lanolin DAB	5651400	1
Henry Lamotte	Lanolin DAB	5860800	1
Henry/Fiebig	Lanolin DAB	5860800	3
Henry Lamotte/Noweda	Lanolin DAB	5860800	1
Lamotte	Lanolin DAB	5860800	1
Henry Lamotte/ Noweda	Lanolin DAB	60000746001	1
Henry Lamotte Oils GmbH...	Lanolin DAB	5860800	1
Lamotte / AHD EK 6,75	Lanolin DAB	5860800	1
Henry Lamotte Oils GmbH...	Lanolin DAB	6000746001	2
Lamotte	Lanolin DAB	6000746001	2
Henry Lamotte/ Noweda	Lanolin DAB	6000746001	1
Lamotte GmbH	Lanolin DAB	6000746001	1
AHD 14.03.17 Ek:6,91	Lanolin DAB	6003682006	1
Lamotte	Lanolin DAB	600746001	1
Bombastus	Lanolin DAB	60080101	3
Henry Lamotte/Noweda	Lanolin DAB	6003682006	1
Bombastus	Lanolin DAB	266981	1
Fargon	Lanolin DAB	13b12-t01-001087	1
Caelo	Lanolin DAB	15024101	1
Fargon	Lanolin DAB	14k11-t01-016976	1
Caelo	Lanolin DAB	16054702	1
Sanacorp, 25.2.17, 4,30EUR	Lanolin DAB	16302301	1
Sanacorp,04.04.17, 4,30EUR	Lanolin DAB	17018401	1
Bombastus	Lanolin DAB	294968	1
AHD 14.5.14, Preis 16,69	Lanolin DAB		1
Caelo	Lanolin DAB	2913U-01765	1
Sanacorp	Lanolin DAB	13062501	1
Sanacorp	Lanolin DAB	13227601	2
Caelo	Lanolin DAB	13326701-3455311103	1
Bombastus	Lanolin DAB	277968	1

- 89 763 spectra from 500 *Apo-Ident* customers from a total of 26 715 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Lanolin DAB* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Lanolin DAB* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	280	0	35 264
Type B	0	680	55	93 139
Type C	0	1528	61	89 763

The substance/substance group *Lanolin DAB* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.8571 %)
Type B	100.0000 % (> 99.9761 %)	92.5170 % (> 92.1088 %)
Type C	100.0000 % (> 99.6586 %)	96.1611 % (> 95.9723 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31003	32863	10.07	25.78
31439	31439	0.00	38.10
31746	32863	3.68	27.16
31794	31794	0.00	26.70
31795	31794	1.89	27.50
31798	31794	4.43	28.42
32863	32863	0.00	28.31

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Lauromacrogols
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30765-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Lauromacrogols; Lauromacrogol 400; Macrogol 20 glycerol monostearate type II; Macrogol 20 glyceroli monostearas; Macrogol 4 lauryl ether; Macrogolglycerol hydroxystearate 40; Macrogolglyceroli - hydroxystearas; Macrogoli aether laurilicus; Macrogoli aether laurilicus 9; Polidocanolum; Polyethylenglycolum laurylicum; Polyoxyaethylenglycerinum monostearinicum; Polyoxyaethylenglycerinum trioxystearinicum; Polysorbate 20; Polysorbate 60; Polysorbatum 20; Polysorbatum 60

Special notes

When selecting the *Lauromacrogols* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Lauromacrogols*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Macrogolglycerol...	12L06-N02	30765	32	1405328

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Lauromacrogol 400	13095802	31025	40	20130523
Caelo	Polysorbate 60	12133619	31033	40	20120531
Fagron	Macrogol 20 glyc...	13H09-N06	31182	40	1407195
Caelo	Polysorbate 20	14019413	31323	40	AR-15-FG-008174-01
Caelo	Polysorbate 60	14027102	31324	40	AR-15-FG-008175-01
Caelo	Lauromacrogol 400	13340701	31372	40	AR-14-FG-012474-01
Caelo	Macrogol 4 laury...	14054401	31381	40	1507230
Caelo	Macrogol 20 glyc...	14033701	31469	40	20140214
Caelo	Macrogol 4 laury...	14225301	31583	40	20140820
Caelo	Macrogolglycerol...	14283805	31608	40	20141014
Caelo	Macrogolglycerol...	14283808	31625	40	AR-15-FG-010151-01
Caelo	Polysorbate 20	14365105	31697	40	20150116
Euro OTC	Lauromacrogol 400	1510024	31768	40	20151102*
Euro OTC	Lauromacrogol 400	1605006	32283	40	20160531*
Fagron	Lauromacrogol 400	16E18-B04-326089	32824	40	20160527
Caelo	Polysorbate 60	170258	33058	40	20170313*
Caelo	Macrogol 4 laury...	16100205	33216	40	20160527
Caelo	Macrogol 20 glyc...	162800	33323	40	20161026*
Caelo	Macrogol 20 glyc...	170801	33324	40	20170428*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 792 spectra of 20 reference samples from the substance/substance group *Lauromacrogols*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 20 different batches.
- 34 752 spectra from a total of 820 batches from further 206 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 1648 spectra from 44 reference samples of the substance/substance group *Lauromacrogols*.
- These include spectra of independent samples from 31 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Polysorbate 20	81979119	30050	40
zu erfassen	Macrogol 4 lauryl ether	zu erfassen	30221	40
Caelo	Polysorbate 60	10037208	30230	40
zu erfassen	Polysorbate 20	zu erfassen	30256	40
Caelo	Macrogol 4 lauryl ether	11014502	30278	40
Caelo	Polysorbate 20	1117205	30292	50
Fagron	Macrogolglycerol hydroxystea...	11C18-N08	30343	40
Fagron	Macrogol 20 glycerol monoste...	11E06-N07	30348	30
Caelo	Lauromacrogol 400	10228209	30364	30
zu erfassen	Macrogolglycerol hydroxystea...	zu erfassen	30503	40
Caelo	Macrogol 20 glycerol monoste...	10326003	30526	40
Caelo	Macrogol 4 lauryl ether	12057804	30530	40
Fagron	Lauromacrogol 400	12E03-N04	30573	40
Fagron	Lauromacrogol 400	12E03-N05	30594	40
Caelo	Macrogol 20 glycerol monoste...	14033701	31469 [†]	20
Caelo	Macrogol 4 lauryl ether	14225301	31583 [†]	20
Caelo	Macrogolglycerol hydroxystea...	14283805	31608 [†]	20
Caelo	Macrogolglycerol hydroxystea...	14283808	31625 [†]	20
Fagron	Macrogolglycerol hydroxystea...	14283808	31640	60
Caelo	Polysorbate 20	14365105	31697 [†]	20
Euro OTC	Lauromacrogol 400	1510024	31768 [†]	18
Caelo	Polysorbate 60	161771	32329	40
Caelo	Polysorbate 20	162335	32330	40
Caelo	Macrogol 20 glycerol monoste...	162800	32702	40
Euro OTC	Lauromacrogol 400	1610034	32816	40
Caelo	Polysorbate 60	170999	33322	40
Caelo	Macrogol 20 glycerol monoste...	161883	33325	40
Caelo	Macrogol 20 glycerol monoste...	161884	33326	40
Euro OTC	Lauromacrogol 400	1703033	33448	40
Caelo	Macrogol 4 lauryl ether	161002	33690	40
Caelo	Macrogol 20 glycerol monoste...	161883	33722	40
Caelo	Macrogol 20 glycerol monoste...	161884	33723	40
Euro OTC	Lauromacrogol 400	1710026	33844	40
Caelo	Polysorbate 20	172825	33884	40
Caelo	Polysorbate 60	172920	33888	40
Caelo	Lauromacrogol 400	172509	33899	40
Caelo	Macrogol 20 glycerol monoste...	172882	33901	40
Caelo	Macrogol 20 glycerol monoste...	172880	33902	40
Caelo	Macrogolglycerol hydroxystea...	16349304	33935	40
Caelo	Macrogol 4 lauryl ether	17028303	33950	40
Euro OTC	Lauromacrogol 400	1801043	34007	40
Caelo	Polysorbate 60	180292	34081	40
Caelo	Lauromacrogol 400	181755	34235	40
Caelo	Macrogol 20 glycerol monoste...	181870	34239	40

- 92226 spectra from a total of 2129 batches from further 431 substances. These spectra were

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2733 spectra from 415 *Apo-Ident* customers from 585 batches from the substance/substance group *Lauromacrogols*.
- These include spectra of independent samples from 564 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macroglycerol hydroxystea...	13072004	11
Caelo	Macroglycerol hydroxystea...	14196805	13
Caelo	Macroglycerol hydroxystea...	14283813	11
Caelo	Macroglycerol hydroxystea...	14283808	7
Fagron	Macroglycerol hydroxystea...	14K28-B01	6
Fagron	Macroglycerol hydroxystea...	14K28-B01-301467	12
Caelo	Macroglycerol hydroxystea...	14C28-b01	1
Fagron	Macroglycerol hydroxystea...	14K28B01301467	1
Caelo	Macroglycerol hydroxystea...	14k28-b01-301467	1
Fagron	Macroglycerol hydroxystea...	14K28-B01-301468	8
Caelo	Macroglycerol hydroxystea...	15169404	10
Caelo	Macroglycerol hydroxystea...	15169405	12
Caelo	Macroglycerol hydroxystea...	15169805	1
Alliance H.	Macroglycerol hydroxystea...	15378302	1
Caelo	Macroglycerol hydroxystea...	15378302	16
Caelo	Macroglycerol hydroxystea...	15378304	5
Caelo	Macroglycerol hydroxystea...	16045403	14
Caelo	Macroglycerol hydroxystea...	16247401	17
Calo/AHD	Macroglycerol hydroxystea...	16247401	1
Caelo	Macroglycerol hydroxystea...	16247405	4
Caelo	Macroglycerol hydroxystea...	16247410	3
Caelo	Macroglycerol hydroxystea...	16349303	3
Caelo	Macroglycerol hydroxystea...	1647401	1
Fagron	Macroglycerol hydroxystea...	16H05-B04-324748	4
Caelo	Macroglycerol hydroxystea...	2101415	1
Caelo	Macroglycerol hydroxystea...	12158606	10
Caelo	Macroglycerol hydroxystea...	12158602	2
Caelo	Macroglycerol hydroxystea...	12071402	1
Caelo	Macroglycerol hydroxystea...	13072005	1
Caelo	Macroglycerol hydroxystea...	13072008	4
Sanacorp	Macroglycerol hydroxystea...	14B14B06	1
Fagron	Macroglycerol hydroxystea...	14B14-B06	8
Fagron	Macroglycerol hydroxystea...	14B14-B06-292160	12
Fagron	Macroglycerol hydroxystea...	14K28-801-301467	1
Fagron	Macroglycerol hydroxystea...	14B14-B06-295486	10
Caelo	Macroglycerol hydroxystea...	29091402	1
Caelo	Macroglycerol hydroxystea...	30091404	1
BASF SE/ Phönix	Macroglycerol hydroxystea...	39658724U0	1
Fagron	Macroglycerol hydroxystea...	49426297V0	2
Phönix	Macroglycerol hydroxystea...	11161504	1
Caelo	Macroglycerol hydroxystea...	6021304	1
Caelo	Macroglycerol hydroxystea...	11161504	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogolglycerol hydroxystea...	12103803	5
Kehr	Macrogolglycerol hydroxystea...	11K16-n09	1
Caelo	Macrogolglycerol hydroxystea...	12103807	5
Anzag 07.02.13 ek:4,70EUR	Macrogolglycerol hydroxystea...	12103807	1
Caelo	Macrogolglycerol hydroxystea...	1215614	1
Caelo	Macrogolglycerol hydroxystea...	12158610	2
Caesar &Loretz GmbH	Macrogolglycerol hydroxystea...	12158610	1
Caelo	Macrogolglycerol hydroxystea...	12158611	4
Caelo	Macrogolglycerol hydroxystea...	12158614	4
Fagron	Macrogolglycerol hydroxystea...	12E11-N03	4
Fagron	Macrogolglycerol hydroxystea...	13D23-N03	16
Fagron	Macrogolglycerol hydroxystea...	12L06-N02	14
Fagron	Macrogolglycerol hydroxystea...	13D23N03	3
Großhandel	Macrogolglycerol hydroxystea...	13D23-N03	1
Fagron	Macrogolglycerol hydroxystea...	13K07-N02	7
Caelo	Macrogolglycerol hydroxystea...	17111202	1
Caelo	Macrogolglycerol hydroxystea...	19111202	1
Sanacorp	Macrogolglycerol hydroxystea...	2/AKD16	1
Caelo	Macrogolglycerol hydroxystea...	20130912-1	1
Caelo	Macrogolglycerol hydroxystea...	21121203	1
Caelo	Macrogolglycerol hydroxystea...	26061303	1
Caelo	Macrogolglycerol hydroxystea...	26061304	1
Fagron	Macrogolglycerol hydroxystea...	87712168E0	1
Spinnrad	Macrogol 4 lauryl ether	12124000321	5
Noweda	Macrogol 4 lauryl ether	512400031	1
Spinnrad/Noweda	Macrogol 4 lauryl ether	512400031	1
Gepepharm GmbH.	Lauromacrogol 400	130109	1
Caelo	Lauromacrogol 400	13340701	7
Caelo	Lauromacrogol 400	13095817	4
Caelo	Lauromacrogol 400	13340703	4
Großhandel	Lauromacrogol 400	13340703	1
Caelo	Lauromacrogol 400	13340706	1
Caelo	Lauromacrogol 400	13340710	3
Sanacorp	Lauromacrogol 400	13340710	1
Caelo	Lauromacrogol 400	13340712	1
Fagron	Lauromacrogol 400	13J01-N03	2
Fagron	Lauromacrogol 400	13K18-B30	5
Caelo	Lauromacrogol 400	140911	8
Caelo	Lauromacrogol 400	14220805	2
Caelo	Lauromacrogol 400	14220807	2
Caelo	Lauromacrogol 400	14220806	1
Caelo	Lauromacrogol 400	14220808	1
Caelo	Lauromacrogol 400	14220813	3
Caelo	Lauromacrogol 400	14220816	8
Caelo	Lauromacrogol 400	14220818	1
Fagron	Lauromacrogol 400	14E06-B05-293846	4
Fagron	Lauromacrogol 400	14E06-B05-293847	12
Caelo	Lauromacrogol 400	14220820	1
Fagron	Lauromacrogol 400	14H04-B03-301517	7
Fagron	Lauromacrogol 400	14H04-B03-301519	6
Caelo	Lauromacrogol 400	15192202	1
Caelo	Lauromacrogol 400	15192203	2
Caelo	Lauromacrogol 400	15192205	8
Caelo	Lauromacrogol 400	15192208	5
Lamotte	Lauromacrogol 400	15192209	1
Caelo	Lauromacrogol 400	15192209	2
Caelo	Lauromacrogol 400	15363702	6
Caelo	Lauromacrogol 400	15363701	1
Caelo	Lauromacrogol 400	15363703	3
Caelo	Lauromacrogol 400	15363707	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lauromacrogol 400	15363710	1
Caelo	Lauromacrogol 400	15363712	6
Caesar & Loretz/Sanacorp	Lauromacrogol 400	15363712	1
Devesa/NOW	Lauromacrogol 400	15363712	1
Fagron	Lauromacrogol 400	15A26-B02	8
apomix PKH GmbH Halle	Lauromacrogol 400	15A26B02-30362	1
Fagron	Lauromacrogol 400	15A26-B02-303628	19
Caelo	Lauromacrogol 400	15a26-b02-303628	1
gepepharm	Lauromacrogol 400	160312	15
gepepharm/Phoenix	Lauromacrogol 400	160312	3
Gepepharm/AHD	Lauromacrogol 400	160312	1
Gepeph/Phönix	Lauromacrogol 400	160312	1
Gepepharm	Lauromacrogol 400	160312	5
gepepharm/ Noweda	Lauromacrogol 400	160312	1
Sanacorp	Lauromacrogol 400	160312	1
gepepharm/AHD	Lauromacrogol 400	160312	1
Noweda	Lauromacrogol 400	160312A	1
Euro OTC	Lauromacrogol 400	1605006-02	5
gepepharm	Lauromacrogol 400	160312A	1
Caelo	Lauromacrogol 400	16102902	8
Caelo	Lauromacrogol 400	16102901	1
Caelo	Lauromacrogol 400	16102904	1
Caelo	Lauromacrogol 400	16102908	1
Caelo	Lauromacrogol 400	16102911	1
Caelo	Lauromacrogol 400	16102912	1
Caelo	Lauromacrogol 400	16102914	2
Caelo	Lauromacrogol 400	16102920	2
Fagron	Lauromacrogol 400	16E-18-B04-321109	1
Fagron	Lauromacrogol 400	16E18-B04	8
Fagron	Lauromacrogol 400	16E18-B04-322722	7
Fagron	Lauromacrogol 400	16E18-B04-324095	8
Fagron	Lauromacrogol 400	16e18-b04-324095	1
Caelo	Lauromacrogol 400	17046401	4
Caelo	Lauromacrogol 400	17046404	1
Caelo	Lauromacrogol 400	17046406	1
Fagron	Lauromacrogol 400	16E18-B04-326089	5
Caelo	Lauromacrogol 400	17046408	3
Caelo	Lauromacrogol 400	17046411	1
Noweda	Lauromacrogol 400	17-0845	1
Caelo	Lauromacrogol 400	17250902	2
Caelo	Lauromacrogol 400	17250904	1
Fagron	Lauromacrogol 400	17A23-B01-335355	3
Fagron	Lauromacrogol 400	17H21-B07-341795	3
Fagron	Lauromacrogol 400	18-0103	1
Caelo	Lauromacrogol 400	11107901	1
Caelo	Lauromacrogol 400	B14-1414	1
Caelo	Lauromacrogol 400	12092501	2
Caelo	Lauromacrogol 400	12092512	3
Phönix Gotha	Lauromacrogol 400	12092516	1
Caelo	Lauromacrogol 400	12092516	4
Caelo	Lauromacrogol 400	12092525	3
Caelo	Lauromacrogol 400	12092530	3
Fagron	Lauromacrogol 400	12A26-N08	9
Fagron	Lauromacrogol 400	12E03-N04	24
gepepharm/Phoenix	Lauromacrogol 400	130109	4
GEHE	Lauromacrogol 400	130109	1
Caelo	Lauromacrogol 400	13095801	3
Caelo	Lauromacrogol 400	13095802	17
Sanacorp	Lauromacrogol 400	13095802	1
Caelo	Lauromacrogol 400	13095809	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lauromacrogol 400	13095811	5
Fagron	Lauromacrogol 400	13A31-N04	3
Anzag	Lauromacrogol 400	13A31N04	1
Fagron	Lauromacrogol 400	13G02-N03	15
Caelo	Macrogol 20 glycerol monoste...	13213601	1
Caelo	Macrogol 20 glycerol monoste...	13213602	3
Caelo	Macrogol 20 glycerol monoste...	14033701	2
Caelo	Macrogol 20 glycerol monoste...	14302001	2
Euro OTC	Macrogol 20 glycerol monoste...	15410502	1
Caelo	Macrogol 20 glycerol monoste...	15410502	1
Caelo	Macrogol 20 glycerol monoste...	15410504	2
Caelo	Macrogol 20 glycerol monoste...	16188401	2
Caelo	Macrogol 20 glycerol monoste...	17080101	2
Caelo	Macrogol 20 glycerol monoste...	18018506	1
Bombastus	Macrogol 20 glycerol monoste...	700906	1
Fagron	Macrogol 20 glycerol monoste...	12H23-N03	4
Caelo	Macrogol 20 glycerol monoste...	12326201	2
Fagron	Macrogol 20 glycerol monoste...	13H09-N06	5
Bombastus	Macrogol 20 glycerol monoste...	13H09-N06	1
Fagron	Macrogol 20 glycerol monoste...	ES61A85388	1
Caelo	Macrogol 20 glycerol monoste...	ES62A93607	1
Fagron	Macrogolglycerol hydroxystea...	267209T0	1
Fagron	Macrogolglycerol hydroxystea...	13K07N02	3
Caelo	Macrogolglycerol hydroxystea...	140415C	1
Gehe	Macrogolglycerol hydroxystea...	14196804	1
phönix	Macrogolglycerol hydroxystea...	14196805	1
Sanacorp	Macrogolglycerol hydroxystea...	14196805	1
Bombastus	Macrogolglycerol hydroxystea...	14196805	1
Caelo	Macrogolglycerol hydroxystea...	14283805	6
Sanacorp	Macrogolglycerol hydroxystea...	1423805	1
Alliance health	Macrogolglycerol hydroxystea...	14283808	1
Caelo	Macrogolglycerol hydroxystea...	14283810	1
Phönix	Macrogolglycerol hydroxystea...	14283813	1
Euro OTC	Macrogolglycerol hydroxystea...	1483808	1
Caelo	Macrogolglycerol hydroxystea...	14B14-B06-292160	1
AHD	Macrogolglycerol hydroxystea...	14B14-B06-295486	1
Fagron	Macrogolglycerol hydroxystea...	14B14-B06-2965486	1
Fagron	Macrogolglycerol hydroxystea...	14K28-B01-309197	2
Caelo	Macrogol 4 lauryl ether	15041704	7
Caelo	Macrogolglycerol hydroxystea...	14K28-B01-309197	1
Caelo	Macrogol 4 lauryl ether	15041801	5
Caelo	Macrogolglycerol hydroxystea...	14K28-B01-301467	2
Bombastus	Macrogolglycerol hydroxystea...	14K28-B01-309197	1
Caelo	Macrogol 4 lauryl ether	15323701	16
Euro OTC	Macrogol 4 lauryl ether	15323704	1
phönix	Macrogolglycerol hydroxystea...	15378302	1
Sanacorp	Macrogolglycerol hydroxystea...	15378309	1
Caelo	Macrogolglycerol hydroxystea...	15378309	2
Fagron	Macrogolglycerol hydroxystea...	15E18-B10	3
Fagron	Macrogolglycerol hydroxystea...	15F18B10-309198	1
Sanacorp/Fischar	Macrogolglycerol hydroxystea...	16045403	1
Euro OTC	Macrogol 4 lauryl ether	1605006-02	1
Caelo	Macrogol 4 lauryl ether	16100201	2
Caelo	Macrogol 4 lauryl ether	16100205	7
Euro OTC	Macrogol 4 lauryl ether	1610034-01	1
Caelo	Macrogol 4 lauryl ether	16102911	1
Caelo	Macrogolglycerol hydroxystea...	1624401	1
Caesar & Loretz GmbH Hi...	Macrogolglycerol hydroxystea...	16247401	1
Caelo	Macrogolglycerol hydroxystea...	16349304	4
Caelo	Macrogolglycerol hydroxystea...	16349305	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Macrogolglycerol hydroxystea...	16H05-B04	2
Fagron	Macrogolglycerol hydroxystea...	16H05-B04-328700	4
Fagron	Macrogolglycerol hydroxystea...	16H-05-B04-328700	1
GEHE	Macrogol 4 lauryl ether	17028301	2
Caelo	Macrogol 4 lauryl ether	17028303	1
Fagron	Macrogol 4 lauryl ether	170303301	1
Caelo	Macrogol 4 lauryl ether	1710026	1
Fagron	Macrogol 4 lauryl ether	17A23-B01-332867	1
Caelo	Macrogol 4 lauryl ether	1710026-01	1
Fagron	Macrogol 4 lauryl ether	17A23-B01-332868	1
Fagron	Macrogol 4 lauryl ether	17H21-Bo7-341795	1
Caelo	Macrogol 4 lauryl ether	18040101	3
Fagron	Macrogolglycerol hydroxystea...	30378697V0	1
Caelo	Macrogolglycerol hydroxystea...	30378697V0	1
Caelo	Macrogol 4 lauryl ether	4214A-02868	1
Spinnrad	Macrogol 4 lauryl ether	316400031	3
Fagron	Macrogolglycerol hydroxystea...	150056P0	1
Fagron	Macrogolglycerol hydroxystea...	11K16-N09	3
Fagron	Macrogolglycerol hydroxystea...	12D23-N03	1
Fagron	Macrogolglycerol hydroxystea...	12H24-N06	4
Noweda	Macrogolglycerol hydroxystea...	12H24-N06	1
Sanacorp	Macrogolglycerol hydroxystea...	13072004	1
Fagron	Macrogolglycerol hydroxystea...	13D23-N03	2
Fagron	Macrogolglycerol hydroxystea...	13G03-N18	4
Caelo	Macrogol 4 lauryl ether	7121201	1
Caelo	Macrogol 4 lauryl ether	12057804	4
Caelo	Macrogol 4 lauryl ether	12057806	4
Euro OTC	Macrogolglycerol hydroxystea...	28227168E0	1
Caelo	Macrogol 4 lauryl ether	12180202	4
Caelo	Macrogol 4 lauryl ether	12180204	3
Anzag	Macrogol 4 lauryl ether	12180204	1
Caelo	Macrogol 4 lauryl ether	13041706	3
Fagron	Macrogol 4 lauryl ether	13136202	1
Fagron	Macrogolglycerol hydroxystea...	11/26-N11	1
Caelo	Macrogolglycerol hydroxystea...	21071201	1
Fagron	Macrogolglycerol hydroxystea...	14B14-B06292160	2
BASF SE	Macrogolglycerol hydroxystea...	14B14B06292160	1
Caelo	Macrogolglycerol hydroxystea...	14b14-b06-2921600855	1
Fagron	Macrogolglycerol hydroxystea...	14B14-B06-295468	1
Caelo	Macrogolglycerol hydroxystea...	14K28-B01	1
Fagron/VDL	Macrogolglycerol hydroxystea...	14K28-B01-301467	1
Caelo	Macrogolglycerol hydroxystea...	14K28-B01-301468	1
Fagron	Macrogolglycerol hydroxystea...	15E18-B10-309198	2
Gehe	Macrogolglycerol hydroxystea...	16247401	1
Sanacorp	Macrogolglycerol hydroxystea...	16247405	1
Caelo	Macrogolglycerol hydroxystea...	16349310	1
Fagron	Macrogolglycerol hydroxystea...	16H05-B04-326362	5
Fagron	Macrogolglycerol hydroxystea...	16H05B04328700	1
Fagron	Macrogolglycerol hydroxystea...	170703MS/17D03B04333	1
Caelo	Macrogolglycerol hydroxystea...	17D03-B04-333357	1
Fagron	Macrogolglycerol hydroxystea...	17D03-B04-333357	1
Fagron	Macrogolglycerol hydroxystea...	17D03-B-04-333357	1
Fagron	Macrogolglycerol hydroxystea...	17D03-B04-334770	2
Fagron	Macrogolglycerol hydroxystea...	30378697V0	1
Fagron	Macrogolglycerol hydroxystea...	30378697Vo	1
Caelo	Macrogolglycerol hydroxystea...	4613M-01308	1
Fagron/Anzag	Macrogolglycerol hydroxystea...	49426297V0	1
Fagron	Macrogolglycerol hydroxystea...	70201688Q0	1
Fagron	Macrogolglycerol hydroxystea...	509A372	1
Post-Apotheke	Macrogolglycerol hydroxystea...	11K16-N09	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Macrogolglycerol hydroxystea...	12L06-N02	2
Noweda	Macrogolglycerol hydroxystea...	13d23-N03	1
Fagron	Macrogolglycerol hydroxystea...	13d23-n03	1
Fargon	Macrogolglycerol hydroxystea...	13G03-N18	1
Phönix 31.01.2014	Macrogolglycerol hydroxystea...	13D23N03	1
Gehe	Macrogolglycerol hydroxystea...	13K07-N02	1
Caelo	Macrogolglycerol hydroxystea...	13K07-N02	1
Phönix Gotha	Macrogolglycerol hydroxystea...	282268E0	1
Fagron	Macrogolglycerol hydroxystea...	87712168E0	1
Fagron	Macrogolglycerol hydroxystea...	700362-0001	1
Fagron	Macrogolglycerol hydroxystea...	8771216E0	1
Caelo	Macrogolglycerol hydroxystea...	10236606	1
Noweda	Macrogolglycerol hydroxystea...	11C18-N08	1
Gehe	Macrogolglycerol hydroxystea...	11K16-N09	1
Fagron	Macrogolglycerol hydroxystea...	79793775L0	1
Gehe	Macrogolglycerol hydroxystea...	79793775L0	1
AHD;2,32EUR;19.06.15	Macrogol 4 lauryl ether	1014400032	1
gepepharm	Lauromacrogol 400	13-0108	2
Gepepharm	Lauromacrogol 400	130109	16
gebepfarm N PZN 01915345	Lauromacrogol 400	130109	1
gepepharm GmbH	Lauromacrogol 400	130109	29
gepepharm	Lauromacrogol 400	130109	50
Fagron	Lauromacrogol 400	13k18-b30-288163	1
Fagron	Lauromacrogol 400	13K18-B30-288165	2
Fagron	Lauromacrogol 400	14H04-B03	7
Fagron	Lauromacrogol 400	14H04-B03-297073	3
VDL;7,75EUR; 30.09.2015	Macrogol 4 lauryl ether	15041801	1
Sanacorp, 21.03.17, 9,45EUR	Macrogol 4 lauryl ether	16100205	1
Fagron	Lauromacrogol 400	13F13-N07	5
Phönix	Lauromacrogol 400	13K18-B30-28816	2
Euro OTC	Lauromacrogol 400	1142238	2
GEHE	Macrogol 4 lauryl ether	21021172/5	1
Löwen City-Apotheke Bar...	Lauromacrogol 400	2042015	1
Euro OTC	Lauromacrogol 400	1210035-01	2
Euro OTC	Lauromacrogol 400	1275280	1
Caelo	Lauromacrogol 400	1184359	1
Spinnrad	Macrogol 4 lauryl ether	1213400033	2
gepepharm/Anzag	Lauromacrogol 400	130109	1
Fagron	Lauromacrogol 400	13-0599	2
Euro OTC	Lauromacrogol 400	1306004-01	7
Caelo	Lauromacrogol 400	1306004-01	2
Euro OTC	Lauromacrogol 400	1306004-02	3
Caelo	Macrogol 4 lauryl ether	13229805	6
Fagron	Lauromacrogol 400	13B13-N01	25
Fagron	Lauromacrogol 400	13G02-N030021/1084	1
Fagron	Lauromacrogol 400	13G02N03	1
Fagron	Lauromacrogol 400	13g02-n03	1
Fagron	Lauromacrogol 400	13G02-N030047/1084	1
Fagron	Lauromacrogol 400	13G02-N030515/1084	1
Caelo	Lauromacrogol 400	13K18-B30	1
Fagron	Lauromacrogol 400	13K18-B30288163	1
Fagron	Lauromacrogol 400	13K18-B30-288163	10
Fagron	Macrogol 4 lauryl ether	13K18-B30-288163	1
Caelo	Lauromacrogol 400	13K18-B30-288163	1
Fagron	Lauromacrogol 400	13K18-B30-288164	4
Caelo	Macrogol 4 lauryl ether	14054401	7
Euro OTC	Macrogol 4 lauryl ether	14054401	1
Euro OTC	Lauromacrogol 400	1407007-01	4
Euro OTC	Macrogol 4 lauryl ether	1407007-01	1
Euro OTC	Lauromacrogol 400	1407007-07	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lauromacrogol 400	14-0731	5
Caelo	Macrogol 4 lauryl ether	14225306	3
Euro OTC	Lauromacrogol 400	1410010-01	1
Caelo	Lauromacrogol 400	1410010-10	1
Fagron	Lauromacrogol 400	14371706	1
Fagron	Lauromacrogol 400	14A23B01332867	2
Fagron	Macrogol 4 lauryl ether	14E06-B05-293846	1
Fagron	Lauromacrogol 400	14E06-B05-293848	3
Fagron	Lauromacrogol 400	14E06-B05293848	1
Fagron	Macrogol 4 lauryl ether	14E06-B05-293848	1
Fagron	Lauromacrogol 400	14E06-B05-295979	5
Fagron	Lauromacrogol 400	14E06-B05-29847	1
GEHE 28.11.2014	Lauromacrogol 400	14H04-B03-297070	1
Fagron	Lauromacrogol 400	14H04-B03-297070	5
Fagron	Lauromacrogol 400	14H04-B03-297071	2
Caelo	Lauromacrogol 400	14H04-B03-297071	2
Fagron	Lauromacrogol 400	14H04-B03-297074	1
Fagron	Lauromacrogol 400	14H04-B03-301518	2
Fagron	Lauromacrogol 400	14H-04-B03-301518	2
Fagron	Lauromacrogol 400	14H04-B03-301520	4
Euro OTC	Lauromacrogol 400	150303201	1
Euro OTC	Lauromacrogol 400	1503032-01	3
Euro OTC	Lauromacrogol 400	1503032-02	3
Euro OTC	Macrogol 4 lauryl ether	1503032-02	1
Fagron	Lauromacrogol 400	15-0596	1
Euro OTC	Lauromacrogol 400	1510024-01	5
Euro OTC	Lauromacrogol 400	1510024-02	6
Fagron	Lauromacrogol 400	1510024-02	2
Caelo	Lauromacrogol 400	1510024-02	1
Euro OTC	Lauromacrogol 400	1510025-01	1
Schärer & Schläpfer AG/...	Lauromacrogol 400	15A26-02	1
Caelo	Lauromacrogol 400	15A26-B02	2
Fagron	Macrogol 4 lauryl ether	15A26B02	1
Schärer & Schläpfer AG/...	Lauromacrogol 400	15A26-B02	1
Fagron	Lauromacrogol 400	15a26-B02	1
Fagron	Lauromacrogol 400	15A26-B02303626	1
Euro OTC	Lauromacrogol 400	15A26-B02303628	1
Noweda	Lauromacrogol 400	15A26-B02-303625	1
Fagron	Lauromacrogol 400	15A26-B02-303625	8
Fagron	Lauromacrogol 400	15A26-B02-303626	5
Caelo	Lauromacrogol 400	15A26-B02-303627	1
Gehe/Fragon	Lauromacrogol 400	15A26B02303627	1
Fagron	Lauromacrogol 400	15A26-B02-303627	4
Fagron	Lauromacrogol 400	15A26-B02-303629	4
Fagron	Lauromacrogol 400	15G10-B025-318670	1
Fagron	Lauromacrogol 400	15G10-B05	5
Caelo	Lauromacrogol 400	15G10-B05	1
Caelo	Lauromacrogol 400	15G10-B05-313605	1
Fagron	Lauromacrogol 400	15G10-B05-313605	4
Fagron	Lauromacrogol 400	15G10-B05-314366	3
phönix	Lauromacrogol 400	15g10b05314366	1
Fagron	Lauromacrogol 400	15G10-B05-31457	1
Fagron	Lauromacrogol 400	15G10-B05-315343	2
Fagron	Lauromacrogol 400	15G10-B05-317456	5
Gehe/Fragon	Lauromacrogol 400	15G10B05317456	1
Fagron	Lauromacrogol 400	15G10-B05-31745	1
Fagron	Lauromacrogol 400	15G10-B05317456	2
Fagron	Lauromacrogol 400	15g10-b05-317456	1
Fagron	Lauromacrogol 400	15G10-B05-317457	4
Fagron	Lauromacrogol 400	15g10-b05-317457	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lauromacrogol 400	15G10-B05-318670	2
Fagon	Lauromacrogol 400	15G10-B05-318670	1
Caelo	Lauromacrogol 400	16/077/M	1
Fagron	Lauromacrogol 400	16/149/P	1
Fagron	Lauromacrogol 400	15G1-B05-313605	1
Fagron	Lauromacrogol 400	16-0290	2
fagr	Lauromacrogol 400	16-0290	1
Euro OTC	Lauromacrogol 400	1605006-01	6
Fagron	Lauromacrogol 400	1605006-01	1
Euro OTC	Lauromacrogol 400	160500602	1
Euro OTC	Lauromacrogol 400	1610034-01	4
Fagron	Lauromacrogol 400	1618B04-324095	1
Fagron	Lauromacrogol 400	16E18804-322722	1
Caelo	Lauromacrogol 400	16D18-B04	1
Fagron	Lauromacrogol 400	16E18-804-321110	1
Schärer & Schläpfer AG/...	Lauromacrogol 400	16E18-B04	1
Fagron	Lauromacrogol 400	16E18-B04326089	2
Fagron	Lauromacrogol 400	16e18-b04-321109	1
Fagron	Lauromacrogol 400	16E18-B04328699	1
Fagron	Lauromacrogol 400	16E18-B04-321109	6
Fagron	Lauromacrogol 400	16e18-b04-321110	1
Fagron	Lauromacrogol 400	16E18-B04-322723	5
Caelo	Lauromacrogol 400	16E18-B04-322723	2
Fagron	Lauromacrogol 400	16E18-B04-322724	2
Caelo	Lauromacrogol 400	16E18-B04-324095	2
Fagron	Lauromacrogol 400	16E18-B04-32495	1
Fagron	Lauromacrogol 400	16E18B04326089	1
Caelo	Lauromacrogol 400	16E18B04326089	1
Noweda	Lauromacrogol 400	16E18-B04-326089	1
Fagron	Lauromacrogol 400	16E18-B04-326090	2
Caelo	Lauromacrogol 400	16E18-B04-327811	1
Fargon GmbH & Co.KG	Lauromacrogol 400	16E18-B04-326091	1
Fagron	Lauromacrogol 400	16E18-B04-326091	3
Fagron	Lauromacrogol 400	16E-18-B04-327811	1
Fagron	Lauromacrogol 400	16E18-B04-327811	4
Bombastus	Lauromacrogol 400	16E18-B04-328582	1
Fagron	Lauromacrogol 400	16E18-B04-328582	2
Caelo	Lauromacrogol 400	16E18-B04-328699	1
Fagron	Lauromacrogol 400	16E18-B04-328699	9
Fagron	Lauromacrogol 400	16E18-B04-324095	1
Euro OTC	Lauromacrogol 400	1703033-01	5
Euro OTC	Lauromacrogol 400	1710026-01/180606P	1
Fagron	Lauromacrogol 400	17A23-B01	4
Fagron	Lauromacrogol 400	17A23-B01-330982	2
Fragon/Noweda	Lauromacrogol 400	17A23-B01-330982	1
Fagron	Lauromacrogol 400	17A23-B01-330983	6
Sanacoro W.E:01.03.18 ...	Lauromacrogol 400	17A23-b01-330983	1
Fagron	Lauromacrogol 400	17A23B01-330982	1
Gehe/Fragon	Lauromacrogol 400	17A23B01332868	1
Fagron	Lauromacrogol 400	17A23-B01-332869	3
Frago/Gehe	Lauromacrogol 400	17A23-B01-332869	1
Caelo	Lauromacrogol 400	17A23-B01-33287	1
Fagron	Lauromacrogol 400	17A23-B01-335356	1
Frago/Noweda	Lauromacrogol 400	17H21-B07	2
Fagron	Lauromacrogol 400	17H21-B07-341796	3
Fagron	Lauromacrogol 400	18A18-B98-344799	1
Euro OTC	Lauromacrogol 400	260315/1	1
Spinnrad/Noweda	Macrogol 4 lauryl ether	314400031	1
Fagron	Lauromacrogol 400	330983	1
Fagron	Lauromacrogol 400	9002-92-0	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Lauromacrogol 400	MLE2018625	1
Fagron	Lauromacrogol 400	12/18-N04	5
Fagron	Lauromacrogol 400	11I07-N05	1
Ebert+Jacobi	Lauromacrogol 400	12A26-N08	1
Fagron	Lauromacrogol 400	12-1074	1
Fagron	Lauromacrogol 400	12e03-n04	1
Phönix	Lauromacrogol 400	12E03-N04	1
Fagron	Lauromacrogol 400	12I18-N04	2
Fagron	Lauromacrogol 400	12I18-N04	1
Euro OTC	Lauromacrogol 400	1302014-01	5
gepepharm GmbH/Noweda	Lauromacrogol 400	130109	4
Caelo	Lauromacrogol 400	1302014-01	2
EoRho	Lauromacrogol 400	1302014-01	2
Fagron	Lauromacrogol 400	130599	1
Phönix Gotha	Lauromacrogol 400	13-0599	1
Caelo	Lauromacrogol 400	130600401	1
Fagron	Lauromacrogol 400	13813-N01	3
Caelo	Lauromacrogol 400	13B13-N01	3
Phönix Gotha	Lauromacrogol 400	13B13-N01	1
Noweda	Lauromacrogol 400	13B13-N01	1
Schärer&Schlöpfer AG/Ph...	Lauromacrogol 400	13B13-N01	1
Caesar/Noweda	Lauromacrogol 400	13B13-N01	1
Fagron	Lauromacrogol 400	13B13-No1	1
Fagron	Lauromacrogol 400	13B13-N01	1
Fagron	Lauromacrogol 400	13B-N01	1
Phönix Gotha	Lauromacrogol 400	13G02-N030461/1210	1
Fagron	Lauromacrogol 400	13G02-No3	3
Fagron	Lauromacrogol 400	13K18-B30-28816	1
Fargon	Lauromacrogol 400	13K18-B30-288166	1
Fagron	Lauromacrogol 400	21121201	1
Sanacorp	Lauromacrogol 400		1
Euro OTC	Lauromacrogol 400	221542	1
Caelo	Lauromacrogol 400	3121505	1
Euro OTC	Lauromacrogol 400	502956	1
Euro OTC	Lauromacrogol 400	1009020-02	2
Fagron	Lauromacrogol 400	11D28-N20	1
Fagron	Lauromacrogol 400	11F28-N08	3
Fagron	Lauromacrogol 400	11f28-n08	1
Caelo	Macrogol 4 lauryl ether	616A-02868	1
Caelo	Macrogol 4 lauryl ether	1014400032	3
Noweda 10.07.2015	Macrogol 4 lauryl ether	1014400032	1
Spinnrad	Macrogol 4 lauryl ether	1014400032	1
Sanacorp	Macrogol 4 lauryl ether	1014400032	1
Caelo	Macrogol 4 lauryl ether	14225303	4
Bombastus	Macrogol 4 lauryl ether	14225306	1
Caelo	Macrogol 4 lauryl ether	14225308	3
Caesar&loretz,Phoenix	Macrogol 4 lauryl ether	14225308	1
Caesar & Loretz GmbH	Macrogol 4 lauryl ether	14225308	1
Hedinger	Macrogol 4 lauryl ether	15041704	1
Caelo	Macrogol 4 lauryl ether	15323704	7
GEHE	Macrogol 4 lauryl ether	15323704	1
Caelo	Macrogol 4 lauryl ether	16100207	2
Caelo	Macrogol 4 lauryl ether	17028301	2
AHD/ Spinnrad	Macrogol 4 lauryl ether	314400031	1
Spinnrad / Alliance HC	Macrogol 4 lauryl ether	314400031	1
Sinnrad GmbH,Bad Seegeb...	Macrogol 4 lauryl ether	316400031	1
Fagron	Macrogol 4 lauryl ether	316400031	1
Spinnrad GmbH	Macrogol 4 lauryl ether	316400031	1
Caelo	Macrogol 4 lauryl ether	5041801	1
Spinnrad	Macrogol 4 lauryl ether	417400031	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Spinnrad/Sanacorp	Macrogol 4 lauryl ether	715400031	2
Caelo	Macrogol 4 lauryl ether	715400031	1
spinnrad	Macrogol 4 lauryl ether	715400031	2
Caesar & Loretz/Noweda	Macrogol 4 lauryl ether	715400031	1
Spinnrad/ Noweda	Macrogol 4 lauryl ether	715400031	1
Gehe/Spinnrad	Macrogol 4 lauryl ether	715400031	1
Gehe/ Spinnrad	Macrogol 4 lauryl ether	715400031	1
Spinnrad	Macrogol 4 lauryl ether	715400031	1
Spinnrad/Noweda	Macrogol 4 lauryl ether	715400031	1
Spinnrad/Sanacorp	Macrogol 4 lauryl ether	1014400032	1
Caelo	Macrogol 4 lauryl ether	1013400032	1
Spinnrad GmbH	Macrogol 4 lauryl ether	1013400032	1
Caelo	Macrogol 4 lauryl ether	3272779	1
Spinnrad/Gehe	Macrogol 4 lauryl ether	1014400032	1
Spinnrad / Noweda	Macrogol 4 lauryl ether	1014400032	1
Spinnrad	Macrogol 4 lauryl ether	101440032	1
spinnrad	Macrogol 4 lauryl ether	1213400033	1
Gehe	Macrogol 4 lauryl ether	1213400033	1
Noweda	Macrogol 4 lauryl ether	123400033	1
Caelo	Macrogol 4 lauryl ether	13229808	2
Gehe	Macrogol 4 lauryl ether	14054401	1
Caesar & Loretz	Macrogol 4 lauryl ether	14054401	1
Caelo	Macrogol 4 lauryl ether	14054405	7
Caesar & Loretz GmbH	Macrogol 4 lauryl ether	14054405	1
sanacorp	Macrogol 4 lauryl ether	14054405	1
Caelo	Macrogol 4 lauryl ether	14225301	3
Caesar & Loretz GmbH Hi...	Macrogol 4 lauryl ether	14283001	1
Spinnrad/Sanacorp	Macrogol 4 lauryl ether	314400031	2
Raben-Apotheke Mildenau	Macrogol 4 lauryl ether	314400031	1
Spinnrad GmbH	Macrogol 4 lauryl ether	314400031	1
spinnrad	Macrogol 4 lauryl ether	314400031	1
Soinnrad	Macrogol 4 lauryl ether	11012014B	1
Sanacorp	Macrogol 4 lauryl ether	1013400032	1
Spinnrad	Macrogol 4 lauryl ether	1111400032	2
Spinnrad GmbH / Phönix	Macrogol 4 lauryl ether	413400031	1
Noweda	Macrogol 4 lauryl ether	1111400032	1
Spinnrad GmbH	Macrogol 4 lauryl ether	1111400032	1
Spinnrad	Macrogol 4 lauryl ether	111140032	1
Caesar & Loretz GmbH, 4...	Macrogol 4 lauryl ether	11220306	1
Gehe	Macrogol 4 lauryl ether	11220306	1
Caelo	Macrogol 4 lauryl ether	11220306	1
Spinnrad/Anzag	Macrogol 4 lauryl ether	12124000321	2
Gehe 12.04.2013	Macrogol 4 lauryl ether	12124000321	1
Spinnrad GmbH	Macrogol 4 lauryl ether	12124000321	1
spinnrad	Macrogol 4 lauryl ether	1214000321	1
Gehe	Macrogol 4 lauryl ether	12180204	1
Caelo	Macrogol 4 lauryl ether	13041702	6
Caesar&lorenz,Anzag	Macrogol 4 lauryl ether	13041702	1
AHD, 25.07.2013	Macrogol 4 lauryl ether	13136202	1
Jenne	Macrogol 4 lauryl ether	13136202	1
Caelo	Macrogol 4 lauryl ether	13136202	2
Gehe	Macrogol 4 lauryl ether	13136206	1
Caesar &Loretz GmbH	Macrogol 4 lauryl ether	13229802	1
Caelo	Macrogol 4 lauryl ether	13229801	1
Spinnrad GmbH	Macrogol 4 lauryl ether	2240003	1
Spinnrad	Macrogol 4 lauryl ether	413400031	4
Caelo	Macrogol 4 lauryl ether	20121112-1	1
Spinnrad/Phönix	Macrogol 4 lauryl ether	413400031	1
Spinnrad/Noweda	Macrogol 4 lauryl ether	413400031	1
Noweda	Macrogol 4 lauryl ether	413400031	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Spinnrad gmbh	Macrogol 4 lauryl ether	512400031	1
Spinnrad	Macrogol 4 lauryl ether	512400031	2
spinnrad	Macrogol 4 lauryl ether	512400031	1
Spinnrad / Noweda	Macrogol 4 lauryl ether	512400031	1
Caelo	Lauromacrogol 400	12345	1
Fagron	Lauromacrogol 400	130109	2
Fagron	Lauromacrogol 400	1314A-02384	1
gepepharm GmbH	Lauromacrogol 400	140911	15
Gepepharm/Phönix	Lauromacrogol 400	140911	2
gepepharm	Lauromacrogol 400	140911	18
gepepharm	Lauromacrogol 400	140911	2
Gepepharm GmbH	Lauromacrogol 400	140911	4
gepepharm/Noweda	Lauromacrogol 400	140911	2
gepepharm/G	Lauromacrogol 400	140911	1
Gehe/gepepharm	Lauromacrogol 400	140911	3
Gepepharm	Lauromacrogol 400	140911	15
Gepepharm/Gehe	Lauromacrogol 400	140911	7
gepepharm GmbH / Allian...	Lauromacrogol 400	140911	1
gepepharm GmbH/Phönix	Lauromacrogol 400	140911	1
gepepharm/Sanacorp	Lauromacrogol 400	14-1414	1
Croda Iberica	Lauromacrogol 400	1510024-02	1
Fagron	Lauromacrogol 400	15A26B02303627	1
Fagron	Lauromacrogol 400	15E27-B99-308668	1
Fagron	Lauromacrogol 400	15G10-B05-313606	1
Fagron	Lauromacrogol 400	15G10B05315343	1
Caelo	Lauromacrogol 400	15G10-B05-315343	1
Bombastus	Lauromacrogol 400	15G10-B05-318670	1
Gepepharm/Gehe	Lauromacrogol 400	160312	1
gepepharm GmbH	Lauromacrogol 400	160312	13
GEHE	Lauromacrogol 400	160312	1
Thesit	Lauromacrogol 400	160312	5
Caelo	Lauromacrogol 400	16102903	1
Fagron	Lauromacrogol 400	16E18-B04-321110	2
Fagron	Lauromacrogol 400	16E18-B04-328737	1
Euro OTC	Lauromacrogol 400		1
phönix	Lauromacrogol 400	17a23-b01-330983	1
Euro OTC	Lauromacrogol 400	171116PC/1703033-01	1
Gehe	Lauromacrogol 400	17A23-B01-332867	1
Fagron	Lauromacrogol 400	17A23-B01-332867	1
Fagron	Lauromacrogol 400	17H21-B07	3
Fagron	Lauromacrogol 400	303628	1
Fagron	Lauromacrogol 400	T-190116-1	1
GEPE/Anzag	Lauromacrogol 400	110408	2
Gepepharm /Fiebig	Lauromacrogol 400	labip-g97-110616-01	1
Gepepharm / Holdermann	Lauromacrogol 400	110408	8
Anzag	Lauromacrogol 400	110408	4
Sanacorp	Lauromacrogol 400	110408	3
Gegepharm/N/26.02.13 PZ...	Lauromacrogol 400	110408	1
Caelo	Lauromacrogol 400	110408	14
Gepepharm GmbH	Lauromacrogol 400	110408	7
Phoenix/gepepharm	Lauromacrogol 400	110408	1
Caelo	Lauromacrogol 400	13-0108	1
Gehe	Lauromacrogol 400	130109	1
Caesar & Lorenz	Lauromacrogol 400	13095802	1
gepepharm/Sanac	Lauromacrogol 400	130109	3
Spinnrad GmbH	Macrogol 4 lauryl ether	1014400032	2
Sanacorp	Macrogol 4 lauryl ether	12180204	1
Fagron	Macrogol 4 lauryl ether	13B13N01	1
Gutenberg-Apotheke	Macrogolglycerol hydroxystea...	13072008	1
Hedinger	Macrogolglycerol hydroxystea...	13333702	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Macrogolglycerol hydroxystea...	14B14-B06-29212160	1
Caelo	Macrogolglycerol hydroxystea...	P146/14	1
Fagron	Macrogolglycerol hydroxystea...	153356P0	1
	Macrogolglycerol hydroxystea...	13D23-N03	1
Caelo	Polysorbate 20	1709930	1
Caelo	Polysorbate 20	714E-02525	1
Caelo	Polysorbate 20		2
Caelo	Polysorbate 20	12031605	1
Caelo	Polysorbate 20	13086509	10
Phoenix	Polysorbate 20	13086517	1
Caelo	Polysorbate 20	13086517	14
Caelo	Polysorbate 60	13258902	5
VDL	Polysorbate 20	13297005	1
Caelo	Polysorbate 20	13297005	14
Euro OTC	Polysorbate 20	13297009	1
Caelo	Polysorbate 20	13297009	11
gehe	Polysorbate 20	13297009	1
Caelo	Polysorbate 20	13297010	18
gehe	Polysorbate 20	13297099	1
Fagron	Polysorbate 20	13J18-B02	1
Fagron	Polysorbate 20	13J18-B02-287378	4
Caelo	Polysorbate 20	14019410	7
Sanacorp	Polysorbate 20	14019410	1
Phönix	Polysorbate 20	14019410	1
Caelo	Polysorbate 20	14019413	1
	Polysorbate 20	14019413	1
Caelo	Polysorbate 20	14019418	24
Cealo	Polysorbate 20	14019418	1
Gehe	Polysorbate 60	14027102	1
Caelo	Polysorbate 60	14027102	1
Caelo	Polysorbate 60	14027103	1
Caelo	Polysorbate 20	14041503	1
Caelo	Polysorbate 20	14213211	2
Euro OTC	Polysorbate 20	14213213	1
Caelo	Polysorbate 20	14213213	6
Caelo	Polysorbate 20	14365105	17
VDL;4,39EUR;23.03.15	Polysorbate 20	14365105	1
Fagron	Polysorbate 20	14365105	1
Caelo	Polysorbate 20	14365109	22
Caelo	Polysorbate 20	14365116	24
Phönix	Polysorbate 20	14365116	2
Fagron	Polysorbate 20	14A23-B01	1
Caelo	Polysorbate 20	1465109	1
Caelo	Polysorbate 60	15067401	4
Caelo	Polysorbate 60	15067407	1
Caesar&Loretz GmbH	Polysorbate 60	15067409	1
Caelo	Polysorbate 20	15159004	1
Caelo	Polysorbate 20	1532302	1
Caelo	Polysorbate 20	15323902	15
Caelo	Polysorbate 20	1523902	1
Noweda	Polysorbate 20	15323902	1
Cael Sanacorp 05.02.2016	Polysorbate 20	15323902	1
WE:31.03.2016 Sanacorp	Polysorbate 20	15323902	2
Caesar & Loretz GmbH Hi...	Polysorbate 20	15323902	1
Caelo	Polysorbate 20	15323912	10
Caelo	Polysorbate 20	1532902	2
Caelo	Polysorbate 60	16007302	5
Caelo	Polysorbate 20	16029703	40
Fagron	Polysorbate 20	16029703	1
Euro OTC	Polysorbate 20	16029703	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Polysorbate 20	16114502	3
Fagron	Polysorbate 20	16114502	2
Sanacorp, 16.08.2016, 4...	Polysorbate 20	16114502	1
Caelo	Polysorbate 20	16114503	19
Sanacorp	Polysorbate 20	16114503	1
Ebert-Jacobi	Polysorbate 60	1617106	1
Caelo	Polysorbate 60	16177104	1
Caelo	Polysorbate 20	16233501	21
Caesar & Loretz GmbH	Polysorbate 20	16233501	1
Phönix	Polysorbate 20	16233501	1
Gehe	Polysorbate 20	16233501	1
Caelo	Polysorbate 20	16233506	6
Noweda	Polysorbate 20	16233506	1
Caelo	Polysorbate 20	17021303	9
Caelo	Polysorbate 20	170130P/16233501	1
Caelo	Polysorbate 20	17021310	8
Alliance	Polysorbate 20	17021310	1
AHD	Polysorbate 20	17021310	1
Caelo	Polysorbate 20	17021311	3
Caelo	Polysorbate 20	17325301	2
Caelo	Polysorbate 60	17099903	1
Caelo	Polysorbate 20	17325407	1
Caelo	Polysorbate 20	1912M-02525	1
Fagron	Polysorbate 20	25092014A	1
Caelo	Polysorbate 20	3114I-02525	1
Caelo	Polysorbate 20	5/AKD7	1
Caelo	Polysorbate 20	611poly	1
Caelo	Polysorbate 20	11117201	1
Caelo	Polysorbate 20	11117205	2
Caelo	Polysorbate 20	11179506	2
Caelo	Polysorbate 20	11179511	5
Noweda	Polysorbate 20	11347701	1
Phönix Gotha	Polysorbate 20	11347701	1
Phönix	Polysorbate 20	11347707	1
Caelo	Polysorbate 20	11347702	1
Caelo	Polysorbate 20	11347707	6
VDL	Polysorbate 20	11347707	1
Caelo	Polysorbate 20	11347711	9
Sanacorp	Polysorbate 20	11347711	1
Caelo	Polysorbate 20	11347715	7
Caelo	Polysorbate 20	11347719	10
Fagron	Polysorbate 20	11I15-N13	1
Caelo	Polysorbate 20	12303006	8
Anzag	Polysorbate 20	12303006Caelo	1
Caelo	Polysorbate 20	12303016	14
Noweda	Polysorbate 20	12303016	1
Sanacorp	Polysorbate 20	12345	1
Caelo	Polysorbate 20	13066512	1
Caelo	Polysorbate 20	13086505	3
GEHE	Polysorbate 20	13086505	1
Sanacorp	Polysorbate 20	13086509	2
Caelo	Polysorbate 20	13086512	16
Caelo	Polysorbate 20	13086510	1
Sanacorp	Polysorbate 20	13086517	1
Alliance Health 15.10.2...	Polysorbate 20	13086517	1
Caelo	Polysorbate 20	1329001	1
Caelo	Polysorbate 20	13297001	8
Bombastus	Polysorbate 20	13297001	1
Phönix Gotha	Polysorbate 20	13297005	1
Caelo	Polysorbate 20	19021303	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Polysorbate 20	20210312-1	1
Caelo	Polysorbate 20	20220312-1	1
Caelo	Polysorbate 20	4612A-02525	1
Caelo	Polysorbate 20	2704212525	1
Caelo	Polysorbate 20	2713E2525	1
Caelo	Polysorbate 20	80479168	1
Caelo	Polysorbate 60	12133611	2
Caelo	Polysorbate 60	12133602	1
Caelo	Polysorbate 60	13098006	1
Gehe am 06.02.2014	Polysorbate 60	13258902	1
Caelo	Polysorbate 60	12133612	1
Caelo	Polysorbate 20	11347701	1
Sanacorp	Polysorbate 20	12303006	1
Caelo	Polysorbate 20	13297001-3449311132	1
Caelo	Macrogol 20 glycerol monoste...	3042223	1
gepepharm GmbH	Lauromacrogol 400	1915345	1
gepepharm GmbH	Lauromacrogol 400	8062015D	1
Gepepharm über Phoenix	Lauromacrogol 400	130109	7
Noweda	Lauromacrogol 400	1915345	1
Caelo	Lauromacrogol 400		2
gepepharm/Otto	Lauromacrogol 400	130109	25
gepepharm/Sanacorp	Lauromacrogol 400	130109	7
gepepharm/ Sanacorp 22...	Lauromacrogol 400	130109	1
Allianc. Healthc./ Gepe	Lauromacrogol 400	130109	4
gepepharm/AHD	Lauromacrogol 400	130109	1
Caelo	Lauromacrogol 400	130109	10
gepepharm /Noweda	Lauromacrogol 400	130109	3
Gepepharm / Noweda	Lauromacrogol 400	140911	1
gepepharm GmbH/ NOW	Lauromacrogol 400	140911	1
Gepepharm über Phoenix	Lauromacrogol 400	140911	3
Phoenix	Lauromacrogol 400	140911	2
gepepharm/Phönix	Lauromacrogol 400	140911	5
Gepepharm/Sanacorp	Lauromacrogol 400	140911	9
Gepepharm/Noweda	Lauromacrogol 400	140911	4
Noweda/gepepharm	Lauromacrogol 400	140911	1
Allianc. Healthc./ Gepe	Lauromacrogol 400	140911	2
Depepharm/Sanacorp	Lauromacrogol 400	140911	2
gepepharm/Gehe	Lauromacrogol 400	140911	3
Anzag	Lauromacrogol 400	140911	1
Phönix Gotha	Lauromacrogol 400	140911	1
gepep/Phoenix	Lauromacrogol 400	140911	1
gepepharm/Gehe	Lauromacrogol 400	140911	1
phoenix/noveda	Lauromacrogol 400	140911	2
apomix AHD	Lauromacrogol 400	140911	11
Fiebig/ gepepharm	Lauromacrogol 400	140911	1
Hedinger/ Sanacorp	Lauromacrogol 400	140911	1
gepepharm / Phönix	Lauromacrogol 400	140911	2
gepepharm / Noweda	Lauromacrogol 400	140911	1
Gepepharm/Anzag	Lauromacrogol 400	140911	1
AHD 21.03.16 EK 36,71	Lauromacrogol 400	140911	1
Kehr/gepepharm	Lauromacrogol 400	140911	1
Gepepharm/Sana	Lauromacrogol 400	140911	1
Cepepharm / Gehe	Lauromacrogol 400	140911	3
Ebert & Jacobi 31.05...	Lauromacrogol 400	140911	2
Audor Pharma	Lauromacrogol 400	140911	2
gepharm GmbH	Lauromacrogol 400	140911	1
Gehe	Lauromacrogol 400	140911	1
gepepharm/Gehe	Lauromacrogol 400	14-1414	2
Caelo	Lauromacrogol 400	14-1414	1
Gepepharm/ Gehe	Lauromacrogol 400	14-1414	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gepepharm/Sanacorp	Lauromacrogol 400	14-1414	1
Gatt-Koller	Lauromacrogol 400	1504/05153915	2
Fagron	Lauromacrogol 400	15A26B02-303627	1
Fagron	Lauromacrogol 400	15C10-B05-314366	1
Gepepharm/ Noweda	Lauromacrogol 400	160132	1
gepepharm	Lauromacrogol 400	16-0290	1
Sanacorp	Lauromacrogol 400	160290	1
gepepharm	Lauromacrogol 400	1603102	1
AHD	Lauromacrogol 400	160312	1
Gepepharm/Noweda	Lauromacrogol 400	160312	3
gepepharm GmbH/Phönix	Lauromacrogol 400	160312	4
Caelo	Lauromacrogol 400	160312	6
Gepepharm/Fiebig	Lauromacrogol 400	160312	1
gepepharm GmbH/Kehr	Lauromacrogol 400	160312	7
gepepharm/ Kehr	Lauromacrogol 400	160312	1
Gepepharm/ Gehe	Lauromacrogol 400	160312	2
Gepepharm/Alliance	Lauromacrogol 400	160312	1
Allianc. Health. / gepe...	Lauromacrogol 400	160312	1
gepepharm/Krieger	Lauromacrogol 400	160312	1
gepepharm gmbh	Lauromacrogol 400	160312	1
gepepharm/Noweda	Lauromacrogol 400	160312	2
Bombastus	Lauromacrogol 400	160312	2
Gepepharm/Anzag	Lauromacrogol 400	160312	3
gepepharm/Gehe	Lauromacrogol 400	160312	3
gepepharm GmbH, Gehe	Lauromacrogol 400	160312	1
Gepepharm GmbH / Gehe	Lauromacrogol 400	160312	1
gepepharm/ Sanacorp	Lauromacrogol 400	160312	1
Noweda	Lauromacrogol 400	160312	1
gepepharm /Phoenix	Lauromacrogol 400	160312	1
Gehe/gepepharm	Lauromacrogol 400	160312	1
Gegepharm	Lauromacrogol 400	160312	1
Gepepharm/Kehr	Lauromacrogol 400	160312A	1
Gegepharm/Phönix	Lauromacrogol 400	160312A	1
gepepharm/Phoenix	Lauromacrogol 400	160312A	1
Caelo	Lauromacrogol 400	16102916	21
Caelo	Lauromacrogol 400	163012	1
Fagron	Lauromacrogol 400	16B18-B04-321109	1
Caelo	Lauromacrogol 400	17046412	1
Fagron	Lauromacrogol 400	17A23-B01-332868	2
Fagron	Lauromacrogol 400	2215I-02384	1
Gepepharm/Kehr	Lauromacrogol 400	51123/2Sie	1
gepepharm GmbH	Lauromacrogol 400	3072014A	1
Gepepharm	Lauromacrogol 400	6111203	1
gebepharm	Lauromacrogol 400	813E-02384	1
gebepharm	Lauromacrogol 400	813Q-02384	1
Hersteller: gepepharm G...	Lauromacrogol 400	90807	1
Caelo	Lauromacrogol 400	90807	1
gepepharm/Phönix	Lauromacrogol 400	1101409	1
Gepepharm	Lauromacrogol 400	110408	6
Ichthyol	Lauromacrogol 400	110408	4
Gepepharm / Gehe	Lauromacrogol 400	110408	1
Gehe	Lauromacrogol 400	110408	6
gegepharm	Lauromacrogol 400	110408	2
gepepharm GmbH	Lauromacrogol 400	110408	7
Noweda	Lauromacrogol 400	110408	4
Gepepharm/VDL	Lauromacrogol 400	110408	3
Phoenix	Lauromacrogol 400	110408	1
Gepepharm Anzag	Lauromacrogol 400	110408	18
Gepepharm/Gehe	Lauromacrogol 400	110408	2
Kehr	Lauromacrogol 400	110408	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
gepepharm,Anzag	Lauromacrogol 400	110408	2
Gepepharm / Sanacorp	Lauromacrogol 400	110408	2
Phönix/Gepepharma	Lauromacrogol 400	110408	6
ANZAG	Lauromacrogol 400	110408	2
gepepharm	Lauromacrogol 400	110408	5
gepepharm/Phönix	Lauromacrogol 400	110408	2
gepepharm, Noweda	Lauromacrogol 400	110408	4
Stumpf/gepepharm	Lauromacrogol 400	110408	1
Phoenix (gepepharm)	Lauromacrogol 400	110408	1
Gehe 12.04.2013	Lauromacrogol 400	110408	1
gepepharm/gehe	Lauromacrogol 400	110408	3
Gebepharma/Anzag	Lauromacrogol 400	110408	1
Geopharm	Lauromacrogol 400	110408	1
gepepharmGmbH	Lauromacrogol 400	110408	1
Gepepharm GmbH / Anzag	Lauromacrogol 400	110408	1
gepepharm / Noweda	Lauromacrogol 400	110408	1
gepepharm [®] GmbH	Lauromacrogol 400	110408	1
Gehe gepepharmn	Lauromacrogol 400	110408	1
Euro OTC	Lauromacrogol 400	1109011-01	2
Caelo	Lauromacrogol 400	11107804	1
Caelo	Lauromacrogol 400	11107905	1
Caelo	Lauromacrogol 400	11107906	1
Fagron	Lauromacrogol 400	12-0719	1
Caelo	Lauromacrogol 400	12092509	1
Fagron	Lauromacrogol 400	1218n04	1
Fagron	Lauromacrogol 400	12A26-N08	1
Caelo	Lauromacrogol 400	12A26-N08	1
gepepharm GmbH	Lauromacrogol 400	13-0108	1
Caelo	Lauromacrogol 400	12E03-N04	1
gepepharm / Sanacorp	Lauromacrogol 400	13-0108	1
Sanacorp	Lauromacrogol 400	130108	1
Phönix	Lauromacrogol 400	13-0108130109	1
Phoenix	Lauromacrogol 400	13-0108	1
gepepharm/Phönix	Lauromacrogol 400	130109	5
Gebepharma/Anzag	Lauromacrogol 400	130109	1
gepepharm/Noweda	Lauromacrogol 400	130109	2
gepepharm GmbH/ Noweda	Lauromacrogol 400	130109	1
Gepepharm / Gehe	Lauromacrogol 400	130109	6
gepepharm / Sanacorp	Lauromacrogol 400	130109	6
Gepepharm/Gehe	Lauromacrogol 400	130109	5
Phönix/Gepepharma	Lauromacrogol 400	130109	4
gepepharm GmbH	Lauromacrogol 400	13-0109	1
gepepharm GmbH	Lauromacrogol 400	130109	1
gepepharm/sanacorp	Lauromacrogol 400	130109	1
Gepepharm/Noweda	Lauromacrogol 400	130109	12
Phönix/gepepharm	Lauromacrogol 400	130109	3
Gepepharm	Lauromacrogol 400	130109	2
gepepharm/ Noweda	Lauromacrogol 400	130109	6
Sanacorp	Lauromacrogol 400	130109	2
Gepepharm/Phönix	Lauromacrogol 400	130109	5
Gepepharm	Lauromacrogol 400	130109	1
gepepharm, Noweda	Lauromacrogol 400	130109	1
gepepharm	Lauromacrogol 400	130109	1
Gepepharm GmbH / Anzag	Lauromacrogol 400	130109	1
Gepepharm/Sanacorp	Lauromacrogol 400	130109	4
gepe/Gehe	Lauromacrogol 400	130109	1
gepepharm/Phönix	Lauromacrogol 400	130109	1
Phönix Gotha	Lauromacrogol 400	130109	3
gehe	Lauromacrogol 400	130109	2
Gepepharm/V.d.L.	Lauromacrogol 400	130109	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gepepharm/Phoenix	Lauromacrogol 400	130109	2
gepepharm / Noweda	Lauromacrogol 400	130109	6
Gepepharm/ Noweda	Lauromacrogol 400	130109	1
Phoenix	Lauromacrogol 400	130109	1
gepepharm/Gehe	Lauromacrogol 400	130109	3
gepepharm gmbh	Lauromacrogol 400	130109	4
gepepharm GmbH / Allian...	Lauromacrogol 400	130109	1
Gepepharm / GEHE	Lauromacrogol 400	130109	1
gepepharma/Phönix	Lauromacrogol 400	130109	1
Gebepharma GmbH	Lauromacrogol 400	130109	2
Noweda/gepepharm	Lauromacrogol 400	130109	1
AHD 09.12.14	Lauromacrogol 400	130109	2
Gepepharm/Fiebig	Lauromacrogol 400	130109	1
gepepharm/Kehr	Lauromacrogol 400	130199	1
gepepharm	Lauromacrogol 400	13044009	1
Gepepharm	Lauromacrogol 400	13051502	1
Caelo	Lauromacrogol 400	13-0599	1
Fagron	Lauromacrogol 400	13-1827	1
Phönix	Lauromacrogol 400	13B13-N01	1
GEHE	Lauromacrogol 400	13B13-N01	1
Euro OTC	Lauromacrogol 400	13F13-N07	1
Hedinger/Gehe	Lauromacrogol 400	13G02-N03	1
Großhandel	Lauromacrogol 400	13J01-N03	1
Calo/AHD	Lauromacrogol 400	14-0731	1
GEPEPHARM	Lauromacrogol 400	140911	5
Gepepharm GmH	Lauromacrogol 400	140911	5
Phoenix/gepepharm	Lauromacrogol 400	140911	3
GEPEPHARM	Lauromacrogol 400	14-1414	1
Fagron	Lauromacrogol 400	14H04B03	1
Gepepharm/Phoenix	Lauromacrogol 400	1614U2384	1
Gepepharm	Lauromacrogol 400	17111403	1
Gepepharm	Lauromacrogol 400	19101201	1
Gepepharm	Lauromacrogol 400	19121203	1
Gepepharm	Lauromacrogol 400	2002312384	1
gebepharma	Lauromacrogol 400	20030412-1	1
Caelo	Lauromacrogol 400	21011501	1
Gepepharm	Lauromacrogol 400	25061502	1
Gepepharm	Lauromacrogol 400	27111202	1
gepepharm GmbH	Lauromacrogol 400	27112014A	1
Gepepharm	Lauromacrogol 400	30031303	1
Fagron	Lauromacrogol 400	3713U-02384	1
Gepepharm/Phoenix	Lauromacrogol 400	3813E02384	1
Gepepharm	Lauromacrogol 400	4111M-02384	1
gepepharm/Otto	Lauromacrogol 400	65793	1
Caeolo/Gehe	Polysorbate 20	16233501	1
Phönix Gotha	Polysorbate 20	14365116	1
Caelo	Polysorbate 20	3478	1
Caelo	Polysorbate 20	2215I-02525	1
Berkel AHK	Polysorbate 20	14019418	1
Fagron	Polysorbate 20	14A23-B01-294167	1
Caelo	Polysorbate 20	2112E-02525	1
Sanacorp	Polysorbate 20	4/AKD1	1

- 88 619 spectra from 501 *Apo-Ident* customers from a total of 26 431 batches from a further 328 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Lauromacrogols* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Lauromacrogols* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	792	0	34 752
Type B	0	1601	47	92 226
Type C	0	2655	78	88 619

The substance/substance group *Lauromacrogols* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.2424 %)
Type B	100.0000 % (> 99.9761 %)	97.1481 % (> 96.9660 %)
Type C	100.0000 % (> 99.6586 %)	97.1460 % (> 97.0362 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30765	30765	0.00	21.87
31025	32283	1.98	23.64
31033	31324	5.37	15.01
31182	31182	0.00	18.19
31323	31323	0.00	24.77
31324	31324	0.00	14.92
31372	31372	0.00	26.33
31381	31381	0.00	23.49
31469	31182	4.61	20.83
31583	31381	7.46	21.50
31608	30765	12.18	26.15
31625	31625	0.00	28.39
31697	31323	10.11	25.17

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31768	31768	0.00	26.57
32283	32283	0.00	25.30
32824	31324	12.72	38.55
33058	33058	0.00	22.45
33216	31381	8.13	25.84
33323	33323	0.00	18.83
33324	33324	0.00	19.30

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Lavender oil / Clary sage oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30861-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Lavender oil / Clary sage oil; Clary sage oil; Lavender oil; Lavendula officinalis oleum; Oleum lavendula officinalis; Oleum salviae sclarea

Special notes

When selecting the *Lavender oil / Clary sage oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Lavender oil / Clary sage oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Lavender oil	12201206	30861	40	20120721
Fagron	Lavender oil	13D24-N01	31223	40	AR-15-FG-009308-01

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Lavender oil	12201213	31237	30	20120721
Caelo	Lavender oil	12334514	31238	30	20121120
Caelo	Clary sage oil	13334802	31413	40	1509218
Caelo	Lavender oil	13252605	31416	40	20130802
Caelo	Lavender oil	13252207	31422	40	20130802
Euro OTC	Lavender oil	1312008-01	31528	40	AR-15-FG-000438-01
Caelo	Clary sage oil	1431102	31668	40	20141031
Euro OTC	Lavender oil	1507031	31721	40	0000010673
Euro OTC	Lavender oil	1510040	31800	40	20151119*
Euro OTC	Lavender oil	1608026	32710	40	20160902*
Fagron	Lavender oil	15I03-B01-314158	32996	40	20150916
Caelo	Clary sage oil	171087	33911	40	20170427*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 540 spectra of 14 reference samples from the substance/substance group *Lavender oil / Clary sage oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 14 different batches.
- 35 004 spectra from a total of 826 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 460 spectra from 14 reference samples of the substance/substance group *Lavender oil / Clary sage oil*.
- These include spectra of independent samples from 9 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Lavender oil	13252207	31422 [†]	20
Euro OTC	Lavender oil	1312008-01	31528 [†]	20
Caelo	Clary sage oil	1431102	31668 [†]	20
Euro OTC	Lavender oil	1507031	31721 [†]	20
Euro OTC	Lavender oil	1510040	31800 [†]	20
Euro OTC	Lavender oil	1601053	31997	40
Euro OTC	Lavender oil	1606004	32458	40
Caelo	Lavender oil	162922	32684	40
Euro OTC	Lavender oil	1610029	32834	40
Euro OTC	Lavender oil	1701005	32986	40
Caelo	Clary sage oil	16269902	33281	40
Euro OTC	Lavender oil	1706003	33597	40
Euro OTC	Lavender oil	1710038	33846	40
Euro OTC	Lavender oil	1808034	34213	40

- 93 414 spectra from a total of 2153 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 776 spectra from 229 *Apo-Ident* customers from 277 batches from the substance/substance group *Lavender oil / Clary sage oil*.
- These include spectra of independent samples from 267 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Lavender oil	13252603	1
VDI;7,53EUR;13.4.15	Lavender oil	13252611	2
Fagron	Lavender oil	13B22-N02	13
Caelo	Lavender oil	14218324	3
Caelo	Lavender oil	13252616	4
Caelo	Lavender oil	13252621	1
Fagron	Lavender oil	14D14-B02	4
Caelo	Lavender oil	15043504	6
Caelo	Lavender oil	15043505	4
Caelo	Lavender oil	15043516	1
Caelo	Lavender oil	15278102	2
Caelo	Lavender oil	273289	1
Fagron	Lavender oil	4A07-B02	1
Fagron	Lavender oil	12B16-N02	4
ACM	Lavender oil	4683Pos2	1
Caelo	Lavender oil	7111403	1
Caelo	Lavender oil	9229520	2
Primavera	Lavender oil	8508143	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Taoasis	Lavender oil	LOG2014-119488	2
Taoasis	Lavender oil	LOG2014-120278	1
Taoasis	Lavender oil	LOG2014-121513	2
Fagron	Lavender oil	11K2900-0002	1
Caelo	Lavender oil	12334513	12
Fagron	Lavender oil	12H15-N01	11
Fagron	Lavender oil	13/11-B03-289423	1
Taoasis	Lavender oil	130201-115678	1
Taoasis	Lavender oil	130201-116395	2
Taoasis / Sanacorp	Lavender oil	130201-115678	1
Taoasis	Lavender oil	130201-116899	6
Taoasis/Noweda	Lavender oil	130201-117994	1
PHOENIX 05.05.2014	Lavender oil	130202301	1
Caelo	Lavender oil	13133703	12
Noweda	Lavender oil	13133703	2
GEHE 17.05.2014	Lavender oil	13133703	1
Caelo	Lavender oil	13133704	3
Caelo	Lavender oil	13133706	21
Caelo	Lavender oil	13133705	1
Caelo	Lavender oil	13133715	11
Caelo	Lavender oil	13202218	1
Caelo	Lavender oil	13252201	5
Noweda	Lavender oil	13252201	1
Caelo	Lavender oil	13252202	7
Caelo	Lavender oil	13252204	8
Caelo	Lavender oil	13252206	3
Caelo	Lavender oil	13252207	4
Caelo	Lavender oil	13252212	1
Caelo	Lavender oil	13252218	1
Caelo	Lavender oil	13252220	7
Phönix	Lavender oil	13252220	2
AHD 15.10.14 EK 5,76	Lavender oil	13252220	1
Caelo	Lavender oil	13252225	4
Caelo	Lavender oil	13252611	2
Fagron	Lavender oil	13D24-0003	1
Fagron	Lavender oil	13D24N01	3
Fagron	Lavender oil	13D24-N01	4
Fagron	Lavender oil	13E30-B01	4
Fagron	Lavender oil	13E30-B01-285543	2
Fagron	Lavender oil	13E30-B01-285544	5
Caelo	Lavender oil	13E30-B01-285545	1
Fagron	Lavender oil	13e30-b01-285544	1
Fagron	Lavender oil	13G11-N04	13
Fagron	Lavender oil	13J11-B03	4
Fagron	Lavender oil	13G11N04	1
Sanacorp	Lavender oil	13G11-N04	1
Fagron	Lavender oil	13J11B03	1
Fagron	Lavender oil	13J11-B03-289423	5
Fagron	Lavender oil	13J11-B03-289424	1
Fagron	Lavender oil	13J11-BC3-289423	1
Caelo	Lavender oil	13J11-B03-289424	1
Caelo	Lavender oil	1414-B02-302035	1
Caelo	Lavender oil	14218301	2
Fagron	Lavender oil	13J11-B03-289425	1
Caelo	Lavender oil	14218311	11
Caelo	Lavender oil	1421831	1
Caelo	Lavender oil	14218312	1
Caelo	Lavender oil	14218314	9
Caelo	Lavender oil	14218321	1
Caelo	Lavender oil	14218322	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lavender oil	14a07-B02	1
Fagron	Lavender oil	14A07-B02	4
Fagron	Lavender oil	14A07B02	1
Frey & Lau GmbH/Phoenix	Lavender oil	14A07-B02-289955	1
Fagron	Lavender oil	14A07-B02-289955	2
Fagron	Lavender oil	14A07-B02-294931	4
Caelo	Lavender oil	14A07-B02-289955	1
Fagron	Lavender oil	14D14B02	1
Fagron	Lavender oil	14d14-b02	1
Fagron	Lavender oil	14D14-B02-302034	3
Fagron	Lavender oil	14D-B02-302034	1
Fagron	Lavender oil	14D14-B02-302035	1
Fagron	Lavender oil	15/03-B01-314158	2
Euro OTC	Lavender oil	1503006-01	2
Caelo	Lavender oil	15043512	1
Caelo	Lavender oil	15043518	3
Caelo	Lavender oil	15043519	3
Caelo	Lavender oil	15043520	2
Caelo	Lavender oil	15043525	13
Euro OTC	Lavender oil	15043525	2
Caelo	Lavender oil	15043526	6
Caelo	Lavender oil	1510040-01	2
Euro OTC	Lavender oil	1507031-01	1
Caelo	Lavender oil	15278101	7
Caelo	Lavender oil	1536007	1
Caelo	Lavender oil	15363001	5
gehe	Lavender oil	15363019	1
Caelo	Lavender oil	15363302	16
AHD	Lavender oil	15363302	1
Caelo	Lavender oil	15363304	5
Caelo	Lavender oil	15363310	10
Fagron	Lavender oil	15b10-b07	1
Euro OTC	Lavender oil	15B10-B07	1
Fagron	Lavender oil	15B10-B07-304078	1
Fagron	Lavender oil	15B10-B07-309880	1
Caelo	Lavender oil	15C26-B02-310775	1
Fagron	Lavender oil	15CV26-B02-309881	1
Fagron	Lavender oil	15I03-B01-893241	1
Fagron	Lavender oil	15C26-B02-310775	1
Caelo	Lavender oil	16023018	1
Caelo	Lavender oil	16292205	9
Caelo	Lavender oil	17041502	1
Caelo	Lavender oil	17051404	1
Caelo	Lavender oil	2189	1
Bombastus	Lavender oil	275898	1
Taoasis	Lavender oil	22015-122102	1
Bombastus	Lavender oil	279919	6
Bombastus	Lavender oil	284486	2
Bombastus	Lavender oil	287254	3
Noweda	Lavender oil	284691	1
Bombastus	Lavender oil	288979	5
Bombastus	Lavender oil	293661	4
Noweda	Lavender oil	293661	1
Bombastus	Lavender oil	295593	3
Bombastus	Lavender oil	295738	3
InfectoPharm/Noweda	Lavender oil	297733	1
Bombastus	Lavender oil	297733	2
Taoasis/Taoasis	Lavender oil	30560-112844	1
Taoasis	Lavender oil	30560-111202	1
Caelo	Lavender oil	31031501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Taoasis	Lavender oil	43703-122404	1
Taoasis	Lavender oil	4651010-118091	1
Taoasis	Lavender oil	4651010118091	1
Fagron	Lavender oil	736340	1
Caelo	Lavender oil	480435519	1
Fagron	Lavender oil	759243	1
phönix	Lavender oil	759243	1
Fagron	Lavender oil	780162	1
Caelo	Lavender oil	780162	1
Fagron	Lavender oil	794323	2
Fagron	Lavender oil	893241	3
Bergland Pharma / Anzag	Lavender oil	K051367	1
Caelo	Lavender oil	213A-01990	1
Caelo	Lavender oil	3091207	1
Caelo	Lavender oil	1021403	1
Caelo	Lavender oil	3011401	1
Fagron	Lavender oil	612211990	2
Caelo	Lavender oil	11053709	5
Caelo	Lavender oil	11221308	1
Noweda	Lavender oil	11221309	1
Noweda	Lavender oil	11221315	2
Caelo	Lavender oil	11221315	1
Caesar & Loretz GmbH, 4. . .	Lavender oil	11355	1
Fagron	Lavender oil	11F09-N07	1
Fagron	Lavender oil	11K17-N06	2
Fagron	Lavender oil	11L22-N16	4
Phönix	Lavender oil	11L22-N16	1
Caelo	Lavender oil	12056802	4
Caesar & Loretz GmbH, 4. . .	Lavender oil	12056804	1
Caelo	Lavender oil	12056804	2
Caelo	Lavender oil	12056807	1
Caelo	Lavender oil	12056809	16
gehe.01.2013	Lavender oil	12056809	1
Caelo	Lavender oil	12056814	5
Phönix, 16.10.2012	Lavender oil	12056814	1
Noweda	Lavender oil	12056815	1
Caelo	Lavender oil	12056815	1
Caelo	Lavender oil	12056816	3
Caesar & Loretz/ Noweda	Lavender oil	12056816	1
Caelo	Lavender oil	12056821	2
Caelo	Lavender oil	12056820	1
Noweda	Lavender oil	12056821	1
Sanacorp	Lavender oil	12056821	1
Caelo	Lavender oil	12056822	11
Caelo	Lavender oil	12081301	2
Caelo	Lavender oil	12201201	1
Caelo	Lavender oil	12201206	2
Caelo	Lavender oil	12201213	2
Caelo	Lavender oil	12201215	1
Caelo	Lavender oil	123344514	1
Caelo	Lavender oil	12334501	6
Caelo	Lavender oil	12334503	3
Caelo	Lavender oil	12334507	4
Großhandel	Lavender oil	12334507	1
Caelo	Lavender oil	12334510	7
Anzag 10.04.13 ek: 5.83EUR	Lavender oil	12334510	1
Caesar & Loretz GmbH D- . . .	Lavender oil	12334510	2
Sanacorp	Lavender oil	12334513	1
Caelo	Lavender oil	12334514	6
Euro OTC	Lavender oil	12334514	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lavender oil	12334516	2
Caelo	Lavender oil	12334516	5
AHC 22.08.13 EK5,83 25ml	Lavender oil	12334516	1
Caelo	Lavender oil	1234501	1
geile	Lavender oil	12334516	1
Caelo	Lavender oil	1234514	2
Caelo	Lavender oil	123456	1
Caelo	Lavender oil	1245	1
Fagron	Lavender oil	12D05-N02	6
Phönix	Lavender oil	12D05-N02	1
Caelo	Lavender oil	12H15-N01	4
Gehe	Lavender oil	12H15-N01	1
Caelo	Lavender oil	12H15-No1	1
Taoasis	Lavender oil	130201-115396	1
Euro OTC	Lavender oil	1302023-01	1
Fagron	Lavender oil	13097R	1
Caelo	Lavender oil	13071302	1
Caesar & Loretz GmbH, 4...	Lavender oil	13133706	1
Caelo	Lavender oil	1313375	1
Gehe	Lavender oil	13252601	1
Caelo	Lavender oil	13252605	1
Fagron	Lavender oil	13B22N02	2
Caelo	Lavender oil	13B22-No2	1
Fagron	Lavender oil	13B22-N02	1
Fagron	Lavender oil	13D24-N01	1
AHD, 25.11.13	Lavender oil	13G11-N04	1
PVL	Lavender oil	1401110	1
primavera	Lavender oil	1480513	1
Fagron	Lavender oil	15/03-B01-317686	1
Fagron	Lavender oil	1503B01317685	1
Caelo	Lavender oil	15081301	1
Euro OTC	Lavender oil	15100-40-01	1
Euro OTC	Lavender oil	1510040-01	2
Sanacorp WE:15.04.16 EK...	Lavender oil	15278101	1
Sanacorp	Lavender oil	15278102	1
Caelo	Lavender oil	15363004	4
Caelo	Lavender oil	15363007	1
Caelo	Lavender oil	15363303	3
Caelo	Lavender oil	15363305	2
Fagron	Lavender oil	15B10-B07-304079	1
Fagron	Lavender oil	15C26-B02	1
Fagron	Lavender oil	15C26-B02-314159	3
Fagron	Lavender oil	15I03-B01-317686	1
Euro OTC	Lavender oil	160105301	1
Fagron	Lavender oil	1601311990	1
Fagron	Lavender oil	15103-B01-317686	1
Caelo	Lavender oil	16023001	7
Caelo	Lavender oil	16023002	5
Caelo	Lavender oil	16023003	4
Caelo	Lavender oil	16023015	1
Caelo	Lavender oil	16023017	5
Caelo	Lavender oil	16023014	1
Sanacorp WE11.04.17 -15...	Lavender oil	16023018	1
Euro OTC	Lavender oil	1610029-01	2
Caelo	Lavender oil	16101203	1
Caelo	Lavender oil	16292206	6
Caelo	Lavender oil	16292207	1
Caelo	Lavender oil	16292211	4
AHD	Lavender oil	16292214	1
Caelo	Lavender oil	16292216	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lavender oil	16A07-B01	2
Fagron	Lavender oil	16A07-B01-321295	2
Fagron	Lavender oil	16A07-B01-321296	2
Fagron	Lavender oil	16A07-B01-32195	1
Caelo	Lavender oil	16A07-B01-32196	1
Caelo	Lavender oil	16H01-B04	1
Fagron	Lavender oil	16H01-B04-329684	1
Fagron	Lavender oil	16H01-B04-329685	3
Fagron	Lavender oil	16H01-B04-335505	1
Euro OTC	Lavender oil	1701005-01	2
Caelo	Lavender oil	1712070L/16292216	1
GEHE	Lavender oil	17132710	1
Caelo	Lavender oil	17132710	1
Caelo	Lavender oil	17132711	1
Sanacorp WE: 09.04.18 E...	Lavender oil	17132710	1
Caelo	Lavender oil	17132712	1
Jenne	Lavender oil	192900-0003	1
Taoasis	Lavender oil	22015-122995	2
Taoasis	Lavender oil	2015-2395	1
Anzag 15.11.2012	Lavender oil	270560	1
Fiebi 29.09.12 ek:8,95	Lavender oil	271465	1
Phönix, 23.10.2012	Lavender oil	271465	1
Bombastus	Lavender oil	272797	5
Hedinger	Lavender oil	273007	1
Bombastus	Lavender oil	273980	2
Bombastus	Lavender oil	273007	1
Bombastus	Lavender oil	273725	1
Anzag 10.04.13 ek: 9,57EUR	Lavender oil	273725	1
Caelo	Lavender oil	276783	1
Bombastus	Lavender oil	276783	2
Bombastus	Lavender oil	281980	1
Bombastus	Lavender oil	297007	3
Bombastus	Lavender oil	300666	2
Caelo	Lavender oil	30081301	1
L10/Bambastus	Lavender oil	301109	1
Bombastus	Lavender oil	303335	2
Bombastus	Lavender oil	301109	1
Bombastus	Lavender oil	305279	1
Caelo	Lavender oil	305279	2
Taoasis GmbH	Lavender oil	30560	1
Taoasis	Lavender oil	30560-110088	1
Taoasis	Lavender oil	30560-108326	1
Taoasis	Lavender oil	30560-111623	2
Taoasis	Lavender oil	30560-113018	1
Taoasis	Lavender oil	30560-113019	1
Bombastus	Lavender oil	308369	2
Bombastus	Lavender oil	311261	3
Bombastus	Lavender oil	311562	1
Caelo	Lavender oil	4110E-01990	1
Fagron	Lavender oil	5013I1990	1
Noweda	Lavender oil		1
Caelo	Lavender oil	559621	1
Fagron	Lavender oil	652607	1
Phönix	Lavender oil	652607	1
Phönix	Lavender oil		1
Anzag 17.08.2012	Lavender oil	680588	1
Fagron	Lavender oil	867858	1
Caelo	Lavender oil	911624	1
Fagron	Lavender oil	911624	1
Bombastus	Lavender oil		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lavender oil	L12334511	1
VDL;15.12.14;20.19EUR	Lavender oil	12201213	1
Caelo	Lavender oil	13252224	1
Caelo	Lavender oil	13252607	1
Caelo	Lavender oil	15363002	1
Caelo	Lavender oil	15363013	1
Euro OTC	Lavender oil	1601053-01	1
Bombastus	Lavender oil	290711	1
Caelo	Lavender oil		1
Gehe	Lavender oil	11221301	1
Fiebig	Lavender oil	266070	1
Bombastus	Lavender oil	274960	1
Bombastus	Lavender oil	276644	1
Euro OTC	Lavender oil	1210031-01	1

- 90 576 spectra from 501 *Apo-Ident* customers from a total of 26 738 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Lavender oil* / *Clary sage oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Lavender oil* / *Clary sage oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	540	0	35 004
Type B	0	460	0	93 414
Type C	0	747	29	90 576

The substance/substance group *Lavender oil* / *Clary sage oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.8889 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.6957 %)
Type C	100.0000 % (> 99.6586 %)	96.2629 % (> 95.8763 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30861	32710	17.32	60.85
31223	31223	0.00	60.30
31237	32710	17.20	67.81
31238	31721	2.55	57.86
31413	31413	0.00	57.58
31416	32710	17.98	62.23
31422	31223	2.04	60.02
31528	31528	0.00	64.16
31668	33911	4.50	60.31
31721	31721	0.00	60.60
31800	31800	0.00	58.77
32710	32710	0.00	62.76
32996	31721	3.72	52.26
33911	33911	0.00	58.29

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Liniment ointments
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30982-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Liniment ointments; Hydrated liniment SR/DAC (NRF p.40); Linimentum aquosum N (NRF S. 39); Linimentum aquosum SR/DAC (NRF S.40); Linimentum nonionicum aquosum DAC; Non-ionic hydrous liniment SR/DAC (NRF P. 39)

Special notes

When selecting the *Liniment ointments* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Liniment ointments*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Leyh Pharma	Hydrated linimen...	124611	30982	80	20121123
Bombastus	Non-ionic hydrou...	282335	31169	40	1405256
apomix	Hydrated linimen...	73A041113	31170	40	1405106

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Hydrated linimen...	14063201	31438	40	20140513
Caelo	Non-ionic hydrou...	14134601	31476	40	20140516
Caelo	Non-ionic hydrou...	15420901	31827	40	20151221
Bombastus	Non-ionic hydrou...	297176	31830	40	20150909
Caelo	Hydrated linimen...	160217	31938	40	20160127*
Caelo	Hydrated linimen...	153244	31939	40	20150922*
Caelo	Non-ionic hydrou...	153267	31943	40	20150922*
Caelo	Non-ionic hydrou...	153342	31945	40	20151008*
Caelo	Non-ionic hydrou...	160216	31947	40	20160127*
Caelo	Non-ionic hydrou...	161138	32090	40	20160415*
apomix	Hydrated linimen...	08B010216	32305	40	20160203
apomix	Hydrated linimen...	32B090516	32727	40	20160509

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 640 spectra of 15 reference samples from the substance/substance group *Liniment ointments*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 15 different batches.
- 34 904 spectra from a total of 825 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 2440 spectra from 63 reference samples of the substance/substance group *Liniment ointments*.
- These include spectra of independent samples from 37 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
apomix	Hydrated liniment SR/DAC (NR...	105A291210	30240	80
Bombastus	Non-ionic hydrous liniment S...	253524	30282	40
Leyh Pharma	Hydrated liniment SR/DAC (NR...	111441	30283	40
Bombastus	Non-ionic hydrous liniment S...	00262833	30332	40
Bombastus	Non-ionic hydrous liniment S...	00262832	30338	40
apomix	Hydrated liniment SR/DAC (NR...	62B010811	30351	30
apomix	Non-ionic hydrous liniment S...	18A011111	30354	30
Caelo	Non-ionic hydrous liniment S...	15420901	31827 [†]	20
Bombastus	Non-ionic hydrous liniment S...	297176	31830 [†]	20
Caelo	Hydrated liniment SR/DAC (NR...	160217	31938 [†]	20
Caelo	Hydrated liniment SR/DAC (NR...	153244	31939 [†]	20
Caelo	Non-ionic hydrous liniment S...	153267	31943 [†]	20
Caelo	Non-ionic hydrous liniment S...	153341	31944	60
Caelo	Non-ionic hydrous liniment S...	153342	31945 [†]	20
Caelo	Non-ionic hydrous liniment S...	153342	31946	60
Caelo	Non-ionic hydrous liniment S...	160216	31947 [†]	20
Caelo	Hydrated liniment SR/DAC (NR...	160217	32089	40
Caelo	Non-ionic hydrous liniment S...	161139	32091	40
Caelo	Non-ionic hydrous liniment S...	161140	32092	40
Caelo	Hydrated liniment SR/DAC (NR...	161824	32395	40
Caelo	Hydrated liniment SR/DAC (NR...	161824	32396	40
Caelo	Hydrated liniment SR/DAC (NR...	161824	32397	40
Caelo	Non-ionic hydrous liniment S...	161767	32398	40
Caelo	Non-ionic hydrous liniment S...	161767	32399	40
Caelo	Non-ionic hydrous liniment S...	161767	32400	40
Caelo	Non-ionic hydrous liniment S...	161768	32401	40
Caelo	Non-ionic hydrous liniment S...	161768	32402	40
Caelo	Non-ionic hydrous liniment S...	161768	32403	40
Caelo	Non-ionic hydrous liniment S...	162035	32404	40
Caelo	Non-ionic hydrous liniment S...	162035	32405	40
Caelo	Non-ionic hydrous liniment S...	161139	32604	40
Caelo	Non-ionic hydrous liniment S...	162181	32605	40
Caelo	Non-ionic hydrous liniment S...	162181	32606	40
Caelo	Non-ionic hydrous liniment S...	162181	32607	40
Caelo	Non-ionic hydrous liniment S...	162181	32608	40
Caelo	Non-ionic hydrous liniment S...	162182	32609	40
Caelo	Non-ionic hydrous liniment S...	162182	32610	40
Caelo	Non-ionic hydrous liniment S...	162183	32611	40
Caelo	Non-ionic hydrous liniment S...	162183	32612	40
Caelo	Non-ionic hydrous liniment S...	162183	32613	40
Caelo	Non-ionic hydrous liniment S...	162183	32614	40
Caelo	Non-ionic hydrous liniment S...	162584	32615	40
Caelo	Non-ionic hydrous liniment S...	162585	32616	40
Caelo	Non-ionic hydrous liniment S...	162585	32617	40
Caelo	Non-ionic hydrous liniment S...	162585	32618	40
Caelo	Non-ionic hydrous liniment S...	163188	32874	40
Caelo	Non-ionic hydrous liniment S...	171173	33400	40
Caelo	Non-ionic hydrous liniment S...	171444	33401	40
Caelo	Non-ionic hydrous liniment S...	171446	33402	40
Caelo	Non-ionic hydrous liniment S...	171445	33403	40
Caelo	Non-ionic hydrous liniment S...	172088	33569	40
Caelo	Non-ionic hydrous liniment S...	172089	33570	40
Caelo	Non-ionic hydrous liniment S...	172090	33571	40
Caelo	Hydrated liniment SR/DAC (NR...	171989	33572	40
Caelo	Non-ionic hydrous liniment S...	172611	33746	40
Caelo	Non-ionic hydrous liniment S...	172612	33747	40
Caelo	Non-ionic hydrous liniment S...	173125	34013	40
Caelo	Hydrated liniment SR/DAC (NR...	173334	34014	40
Caelo	Non-ionic hydrous liniment S...	180507	34146	40
Caelo	Non-ionic hydrous liniment S...	180867	34147	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Non-ionic hydrous liniment S...	181744	34148	40
Caelo	Hydrated liniment SR/DAC (NR...	181848	34158	40
Caelo	Non-ionic hydrous liniment S...	182331	34254	40

- 91 434 spectra from a total of 2123 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1929 spectra from 275 *Apo-Ident* customers from 779 batches from the substance/substance group *Liniment ointments*.
- These include spectra of independent samples from 772 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Apomix	Non-ionic hydrous liniment S...	23A111213	3
apomix	Hydrated liniment SR/DAC (NR...		1
Apomix	Hydrated liniment SR/DAC (NR...	5A210113	1
hochhaus-Apotheke	Non-ionic hydrous liniment S...	2116617027801	1
GEHE	Non-ionic hydrous liniment S...	2A010217	3
Apomix/Noweda	Non-ionic hydrous liniment S...	2A020216	1
Sanacorp	Hydrated liniment SR/DAC (NR...	2A120115	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	2A120115	1
apomix	Non-ionic hydrous liniment S...	2A100215	1
Apomix	Hydrated liniment SR/DAC (NR...	3A070114	2
apomix	Non-ionic hydrous liniment S...	3A100215	2
apomix / AHD	Hydrated liniment SR/DAC (NR...	3B090117	1
Apomix	Hydrated liniment SR/DAC (NR...	3B130115	1
Apomix PKH GmbH	Non-ionic hydrous liniment S...	4A010316	1
Apomix / Gehe	Non-ionic hydrous liniment S...	6A310316	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	5A190115	1
Gehe	Hydrated liniment SR/DAC (NR...	6B260115	1
Apomix / Phönix	Non-ionic hydrous liniment S...	8A100516	3
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	9A010216	3
apomix/Holdermann	Hydrated liniment SR/DAC (NR...	9A020215	1
Pharmachem GmbH & Co KG	Hydrated liniment SR/DAC (NR...	100315	4
PKH Halle/Holdermann	Non-ionic hydrous liniment S...	10A100517	1
apomix	Non-ionic hydrous liniment S...	10A310516	1
Apomix	Non-ionic hydrous liniment S...	10A310516	2
apomix	Hydrated liniment SR/DAC (NR...	10B160215	1
Apomix / Phönix	Non-ionic hydrous liniment S...	12A100615	1
Leyh-Pharma GmbH	Hydrated liniment SR/DAC (NR...	133721	1
apomix / Anzag	Non-ionic hydrous liniment S...	13A040614	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	13A160216	1
Apomix	Hydrated liniment SR/DAC (NR...	13A110214	1
Apomix	Hydrated liniment SR/DAC (NR...	12B110214	1
PKH Halle	Non-ionic hydrous liniment S...	13A200616	1
AHD/23.07.14/Leyh Pharma	Hydrated liniment SR/DAC (NR...	142231	1
Leyh-Pharma GmbH	Hydrated liniment SR/DAC (NR...	142231	1
Caelo	Non-ionic hydrous liniment S...	14276002	8
Caelo	Hydrated liniment SR/DAC (NR...	14350002	1
AHD/Leyh Pharma/25.10.14	Hydrated liniment SR/DAC (NR...	144031	1
apomix	Non-ionic hydrous liniment S...	14A180716	2
Caelo	Hydrated liniment SR/DAC (NR...	15002302	1
AHD /LeyhPharma	Hydrated liniment SR/DAC (NR...	150321	1
Caelo	Non-ionic hydrous liniment S...	15192401	13
Caelo	Hydrated liniment SR/DAC (NR...	15324402	1
Caelo	Non-ionic hydrous liniment S...	15096002	11
Caelo	Non-ionic hydrous liniment S...	15192502	10
Caelo	Hydrated liniment SR/DAC (NR...	15324403	1
Caelo	Non-ionic hydrous liniment S...	15326702	10
Caelo	Non-ionic hydrous liniment S...	15326703	8
Noweda	Hydrated liniment SR/DAC (NR...	15A040315	1
Leyh Pharma	Hydrated liniment SR/DAC (NR...	154111	1
Leyh Pharma / Alliance ...	Hydrated liniment SR/DAC (NR...	153511	1
Apomix / Phönix	Non-ionic hydrous liniment S...	15A150715	1
Pharmachem/Phönix	Hydrated liniment SR/DAC (NR...	160115	1
Leyhs Pharma GmbH	Hydrated liniment SR/DAC (NR...	161121	1
Caelo	Non-ionic hydrous liniment S...	16114002	7
Caelo	Non-ionic hydrous liniment S...	16176701	13
Caelo	Non-ionic hydrous liniment S...	16176802	3
Ichthyol/Phönix	Non-ionic hydrous liniment S...	16113901	2
Caelo	Hydrated liniment SR/DAC (NR...	16182403	1
Caelo	Non-ionic hydrous liniment S...	16218201	11
Caelo	Non-ionic hydrous liniment S...	16218301	9
Caelo	Non-ionic hydrous liniment S...	1625401	1
Caelo	Non-ionic hydrous liniment S...	16318801	17
Leyh Pharma/Sanacorp	Hydrated liniment SR/DAC (NR...	164641	1
Caelo	Non-ionic hydrous liniment S...	1676802	1
Apomix / Phönix	Non-ionic hydrous liniment S...	16A040814	1
Apomix/Kehr	Hydrated liniment SR/DAC (NR...	16B190716	1
Leyhs	Hydrated liniment SR/DAC (NR...	164541	2
Euro OTC	Non-ionic hydrous liniment S...	16C190716	1
Caelo	Hydrated liniment SR/DAC (NR...	17000560	2
Bombastus	Non-ionic hydrous liniment S...	17001409	1
Bombastus	Non-ionic hydrous liniment S...	17001405	1
Bombastus	Hydrated liniment SR/DAC (NR...		2
Caelo	Non-ionic hydrous liniment S...	17027801	4
Caelo	Non-ionic hydrous liniment S...	17027901	3
Caelo	Non-ionic hydrous liniment S...	17064402	9
Caelo	Non-ionic hydrous liniment S...	17175201	1
Caelo	Non-ionic hydrous liniment S...	17175301	2
Caelo	Non-ionic hydrous liniment S...	17261201	2
Apomix	Non-ionic hydrous liniment S...	17A140817	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	17A030314	1
apomix	Hydrated liniment SR/DAC (NR...	17A290216	1
apomix	Hydrated liniment SR/DAC (NR...	17B270317	1
Caelo	Non-ionic hydrous liniment S...	18050701	2
Apomix/Phönix	Non-ionic hydrous liniment S...	18B010916	1
Pharmachem/Phönix	Hydrated liniment SR/DAC (NR...	190416	1
Pharmachem GmbH Pößneck	Hydrated liniment SR/DAC (NR...	190416	1
Bombastus	Non-ionic hydrous liniment S...	192621	1
Apomix PKH GmbH	Non-ionic hydrous liniment S...	20A280915	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix / Phönix	Non-ionic hydrous liniment S...	21A241016	1
PKH/Noweda	Non-ionic hydrous liniment S...	21A250917	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	22A250314	1
Apomix	Non-ionic hydrous liniment S...	21A250917	1
apomix/Phönix	Hydrated liniment SR/DAC (NR...	22A250314	1
Apomix / Phönix	Non-ionic hydrous liniment S...	22A101114	1
Apomix	Hydrated liniment SR/DAC (NR...	22A180417	1
Apomix	Hydrated liniment SR/DAC (NR...	22A250314	1
Apomix/Kehr	Hydrated liniment SR/DAC (NR...	24A070415	1
Apomix / Phönix	Non-ionic hydrous liniment S...	24A271114	1
Apomix/Anzag	Hydrated liniment SR/DAC (NR...	24A290316	1
GEHE/Apomix/29.04.14	Hydrated liniment SR/DAC (NR...	24A310314	1
Apomix/ Phoenix	Non-ionic hydrous liniment S...	25A081215	1
Caelo	Non-ionic hydrous liniment S...	26A281216	1
Bombastus	Hydrated liniment SR/DAC (NR...	26B0590416	1
Raben-Apotheke Mildenau	Hydrated liniment SR/DAC (NR...	25A310314	2
Bombastus	Hydrated liniment SR/DAC (NR...	281644	1
Bombastus	Hydrated liniment SR/DAC (NR...	284234	1
Bombastus	Non-ionic hydrous liniment S...	282945	2
Bombastus	Hydrated liniment SR/DAC (NR...	282885	2
Bombastus	Non-ionic hydrous liniment S...	282947	5
Bombastus	Hydrated liniment SR/DAC (NR...	285737	4
Bombastus	Hydrated liniment SR/DAC (NR...	284340	4
Bombastus	Hydrated liniment SR/DAC (NR...	284339	2
Bombastus	Hydrated liniment SR/DAC (NR...	285738	3
Bombastus	Non-ionic hydrous liniment S...	285793	5
Bombastus	Hydrated liniment SR/DAC (NR...	285739	2
Bombastus	Hydrated liniment SR/DAC (NR...	285792	1
Bombastus	Non-ionic hydrous liniment S...	286417	5
Bombastus	Hydrated liniment SR/DAC (NR...	285838	1
Bombastus	Non-ionic hydrous liniment S...	286418	9
Bombastus	Non-ionic hydrous liniment S...	286419	3
Bombastus	Non-ionic hydrous liniment S...	286421	1
Bombastus	Hydrated liniment SR/DAC (NR...	286935	4
Bombastus	Hydrated liniment SR/DAC (NR...	286936	2
Phönix Gotha	Non-ionic hydrous liniment S...	287316	1
Bombastus	Non-ionic hydrous liniment S...	287318	7
Bombastus	Non-ionic hydrous liniment S...	287321	7
Bombastus	Non-ionic hydrous liniment S...	287323	4
Bombastus	Hydrated liniment SR/DAC (NR...	287410	4
Bombastus	Hydrated liniment SR/DAC (NR...	287409	1
Fagron	Non-ionic hydrous liniment S...	288452	1
Bombastus	Non-ionic hydrous liniment S...	288453	7
Bombastus	Non-ionic hydrous liniment S...	288391	1
Bombastus	Hydrated liniment SR/DAC (NR...	287411	1
Caelo	Non-ionic hydrous liniment S...	288454	1
Bombastus	Non-ionic hydrous liniment S...	288454	6
Bombastus	Hydrated liniment SR/DAC (NR...	288533	1
Bombastus	Hydrated liniment SR/DAC (NR...	288535	6
Bombastus	Hydrated liniment SR/DAC (NR...	288538	2
Gehe	Non-ionic hydrous liniment S...	289208	1
Bombastus	Non-ionic hydrous liniment S...	289864	3
Bombastus	Non-ionic hydrous liniment S...	289209	4
Phönix Gotha	Non-ionic hydrous liniment S...	290610	1
Bombastus	Non-ionic hydrous liniment S...	290121	2
Bombastus	Non-ionic hydrous liniment S...	290613	8
Bombastus	Non-ionic hydrous liniment S...	291222	6
Bombastus	Non-ionic hydrous liniment S...	291223	4
Bombastus	Non-ionic hydrous liniment S...	291860	3
Bombastus	Non-ionic hydrous liniment S...	291859	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrated liniment SR/DAC (NR...	291866	2
Bombastus	Hydrated liniment SR/DAC (NR...	292491	7
Bombastus	Hydrated liniment SR/DAC (NR...	292492	1
Phönix Gotha	Non-ionic hydrous liniment S...	292617	1
Bombastus	Non-ionic hydrous liniment S...	292742	1
Bombastus	Hydrated liniment SR/DAC (NR...	292621	2
Bombastus	Non-ionic hydrous liniment S...	292618	3
Bombastus	Hydrated liniment SR/DAC (NR...	293590	5
Bombastus	Non-ionic hydrous liniment S...	294211	8
Bombastus	Hydrated liniment SR/DAC (NR...	294433	4
Bombastus	Hydrated liniment SR/DAC (NR...	294434	3
Phönix Gotha	Non-ionic hydrous liniment S...	295072	1
Bombastus	Hydrated liniment SR/DAC (NR...	295854	5
Bombastus	Hydrated liniment SR/DAC (NR...	295857	3
Bombastus	Non-ionic hydrous liniment S...	296106	3
Phönix Gotha	Non-ionic hydrous liniment S...	296108	1
Bombastus	Non-ionic hydrous liniment S...	296110	3
Bombastus	Hydrated liniment SR/DAC (NR...	296503	1
Bombastus	Hydrated liniment SR/DAC (NR...	296504	5
Allianc. Health. /Apomix	Hydrated liniment SR/DAC (NR...	297052	1
Bombastus	Hydrated liniment SR/DAC (NR...	296506	3
Bombastus	Non-ionic hydrous liniment S...	29715	1
Bombastus	Non-ionic hydrous liniment S...	297175	3
Caelo	Non-ionic hydrous liniment S...	297171	1
Bombastus	Hydrated liniment SR/DAC (NR...	299986	1
Bombastus	Hydrated liniment SR/DAC (NR...	299987	3
Allianc. Health. /Apomix	Hydrated liniment SR/DAC (NR...	299338	1
Bombastus	Non-ionic hydrous liniment S...	297177	6
Bombastus	Hydrated liniment SR/DAC (NR...	299988	2
PKH Halle / Noweda	Hydrated liniment SR/DAC (NR...	29A270415	1
Bombastus	Non-ionic hydrous liniment S...	300195	4
Bombastus	Non-ionic hydrous liniment S...	300197	5
Bombastus	Non-ionic hydrous liniment S...	300200	7
Phoenix	Non-ionic hydrous liniment S...	300860	1
Bombastus	Non-ionic hydrous liniment S...	300862	5
Bombastus	Non-ionic hydrous liniment S...	300863	4
Bombastus	Non-ionic hydrous liniment S...	300864	6
Ebert u. Jacobi 16.0...	Non-ionic hydrous liniment S...	300864	1
Bombastus	Hydrated liniment SR/DAC (NR...	300965	3
Bombastus	Hydrated liniment SR/DAC (NR...	300968	5
Bombastus	Hydrated liniment SR/DAC (NR...	300966	2
Bombastus	Hydrated liniment SR/DAC (NR...	300967	1
Allianc. Health. /Apomix	Non-ionic hydrous liniment S...	301558	1
Bombastus	Non-ionic hydrous liniment S...	301558	4
Bombastus	Non-ionic hydrous liniment S...	301557	3
Bombastus	Hydrated liniment SR/DAC (NR...	301562	4
Bombastus	Non-ionic hydrous liniment S...	302287	6
Bombastus	Non-ionic hydrous liniment S...	302286	7
Bombastus	Non-ionic hydrous liniment S...	302851	3
Bombastus	Hydrated liniment SR/DAC (NR...	302930	2
Bombastus	Hydrated liniment SR/DAC (NR...	302931	2
Bombastus	Hydrated liniment SR/DAC (NR...	302933	3
Bombastus	Non-ionic hydrous liniment S...	303893	1
Bombastus	Hydrated liniment SR/DAC (NR...	303961	2
Bombastus	Non-ionic hydrous liniment S...	303894	2
Bombastus	Non-ionic hydrous liniment S...	303897	1
Caelo	Non-ionic hydrous liniment S...	303898	1
Bombastus	Non-ionic hydrous liniment S...	303898	3
Bombastus	Hydrated liniment SR/DAC (NR...	303960	1
Bombastus	Hydrated liniment SR/DAC (NR...	303962	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrated liniment SR/DAC (NR...	303963	5
Bombastus	Hydrated liniment SR/DAC (NR...	304842	3
Bombastus	Hydrated liniment SR/DAC (NR...	304843	2
Noweda	Hydrated liniment SR/DAC (NR...	304844	2
Bombastus	Hydrated liniment SR/DAC (NR...	304844	2
Bombastus	Hydrated liniment SR/DAC (NR...	3048843	1
Bombastus	Non-ionic hydrous liniment S...	304933	6
Bombastus	Non-ionic hydrous liniment S...	304932	5
Bombastus	Non-ionic hydrous liniment S...	305594	1
Bombastus	Non-ionic hydrous liniment S...	304934	4
Bombastus	Non-ionic hydrous liniment S...	305595	3
Bombastus	Non-ionic hydrous liniment S...	305596	3
Bombastus	Hydrated liniment SR/DAC (NR...	306059	4
Bombastus	Non-ionic hydrous liniment S...	306361	7
Bombastus	Non-ionic hydrous liniment S...	306362	6
Caelo	Non-ionic hydrous liniment S...	306362	2
Bombastus	Hydrated liniment SR/DAC (NR...	307263	8
Bombastus	Hydrated liniment SR/DAC (NR...	307264	5
Pharmachem / AHD	Hydrated liniment SR/DAC (NR...	307264	1
Bombastus	Non-ionic hydrous liniment S...	307684	3
Bombastus	Non-ionic hydrous liniment S...	308628	9
Bombastus	Hydrated liniment SR/DAC (NR...	307683	1
Bombastus	Hydrated liniment SR/DAC (NR...	307265	2
Bombastus	Non-ionic hydrous liniment S...	308622	1
Bombastus	Non-ionic hydrous liniment S...	309802	3
Bombastus	Non-ionic hydrous liniment S...	308626	3
Bombastus	Hydrated liniment SR/DAC (NR...	308911	1
AMH Niemann GmbH und Co...	Hydrated liniment SR/DAC (NR...	30B040515	1
Bombastus	Non-ionic hydrous liniment S...	310359	1
Bombastus	Non-ionic hydrous liniment S...	310358	2
Bombastus	Hydrated liniment SR/DAC (NR...	310641	1
Bombastus	Hydrated liniment SR/DAC (NR...	310642	2
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	35A080514	1
Noweda PKH Halle	Hydrated liniment SR/DAC (NR...	31A290517	1
Apomix	Hydrated liniment SR/DAC (NR...	34B080514	1
Apomix/Sanacorp	Hydrated liniment SR/DAC (NR...	35A120617	1
Apomix	Hydrated liniment SR/DAC (NR...	31A290517	1
apomix/Holdermann	Hydrated liniment SR/DAC (NR...	35A170516	1
apomix GmbH/Phönix	Hydrated liniment SR/DAC (NR...	36A120514	1
Phönix	Hydrated liniment SR/DAC (NR...	39A190617	1
apomix	Hydrated liniment SR/DAC (NR...	37B010615	1
apomix	Hydrated liniment SR/DAC (NR...	36A120514	1
apomix PKH GmbH Halle	Hydrated liniment SR/DAC (NR...	41A080615	1
apomix	Hydrated liniment SR/DAC (NR...	428140616	1
Gehe/18.08.14	Hydrated liniment SR/DAC (NR...	42A020614	3
Apomix /Noweda	Hydrated liniment SR/DAC (NR...	43B040614	1
PKH Halle	Hydrated liniment SR/DAC (NR...	42A020614	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	41A260514	1
apomix/Gehe	Hydrated liniment SR/DAC (NR...	43A140616	5
Apomix	Hydrated liniment SR/DAC (NR...	45A170614	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	45A170614	1
Apomix	Hydrated liniment SR/DAC (NR...	45A270606	1
apomix	Hydrated liniment SR/DAC (NR...	46B260514	1
apomix	Hydrated liniment SR/DAC (NR...	46A290615	1
Gehe	Hydrated liniment SR/DAC (NR...	46B260514	2
AMH Niemann/Sanacorp	Hydrated liniment SR/DAC (NR...	46B260514	1
Apomix/Phönix	Hydrated liniment SR/DAC (NR...	47B170614	1
Bombastus	Non-ionic hydrous liniment S...	488451	1
Gehe	Hydrated liniment SR/DAC (NR...	50B300614	1
apomix	Hydrated liniment SR/DAC (NR...	53B140716	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix/NOWEDA	Hydrated liniment SR/DAC (NR...	54A030815	1
Apomix /Noweda	Hydrated liniment SR/DAC (NR...	55A250716	1
Apomix/Sanacorp	Hydrated liniment SR/DAC (NR...	55A040815	1
apomix	Hydrated liniment SR/DAC (NR...	55A250716	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	59A050814	1
Apomix/Phönix	Hydrated liniment SR/DAC (NR...	59A050814	1
AHD PKH	Hydrated liniment SR/DAC (NR...	62A021017	1
apomix	Hydrated liniment SR/DAC (NR...	63A120814	1
Apomix	Hydrated liniment SR/DAC (NR...	64A140915	1
apomix	Hydrated liniment SR/DAC (NR...	66A161017	2
Apomix /Noweda	Hydrated liniment SR/DAC (NR...	66b120916	1
PKH Halle/Kehr	Hydrated liniment SR/DAC (NR...		1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	67A010914	1
Gehe 09.01.15 Apomix	Hydrated liniment SR/DAC (NR...	67A010914	1
apomix/Gehe	Hydrated liniment SR/DAC (NR...	68A080914	1
apomix	Hydrated liniment SR/DAC (NR...	68B210915	1
PKH Halle	Hydrated liniment SR/DAC (NR...	69A080914	1
Apomix/Kehr	Hydrated liniment SR/DAC (NR...	69A200916	1
Apomix	Hydrated liniment SR/DAC (NR...	69A220915	1
Apomix / Phönix	Hydrated liniment SR/DAC (NR...	69B071013	2
apomix/Gehe	Hydrated liniment SR/DAC (NR...	72A270916	4
Apomix	Hydrated liniment SR/DAC (NR...	72B041113	1
Apomix /Noweda	Hydrated liniment SR/DAC (NR...	73A101016	1
Gehe 09.01.15/apomix	Hydrated liniment SR/DAC (NR...	74A131014	1
Apomix	Hydrated liniment SR/DAC (NR...	76A271014	1
Noweda	Hydrated liniment SR/DAC (NR...	74A261015	1
Rabenapotheke	Hydrated liniment SR/DAC (NR...	76A271014	1
Gehe/apomix/12.03.14	Hydrated liniment SR/DAC (NR...	79A091213	1
Sanacorp	Non-ionic hydrous liniment S...	8/AKD101	1
apomix/Noweda	Hydrated liniment SR/DAC (NR...	80A211116	2
Apomix/Gehe/12.05.14	Hydrated liniment SR/DAC (NR...	81A161213	1
apomix	Hydrated liniment SR/DAC (NR...	81B071215	1
Apomix	Hydrated liniment SR/DAC (NR...	80A241114	1
apomix	Hydrated liniment SR/DAC (NR...	83B051216	1
apomix / Noweda	Hydrated liniment SR/DAC (NR...	83B051216	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	84B021214	1
Apomix	Hydrated liniment SR/DAC (NR...	85A281215	1
apomix	Hydrated liniment SR/DAC (NR...	88A161214	1
Phönix	Hydrated liniment SR/DAC (NR...		1
Phönix Gotha	Non-ionic hydrous liniment S...	264180	1
Bombastus	Non-ionic hydrous liniment S...	258820	1
Bombastus	Non-ionic hydrous liniment S...	264181	1
Bombastus	Non-ionic hydrous liniment S...	264191	1
Bombastus	Non-ionic hydrous liniment S...	268125	1
Caelo	Non-ionic hydrous liniment S...	269134	1
ANZAG am 23.04.2012	Non-ionic hydrous liniment S...	269132	1
Bombastus	Non-ionic hydrous liniment S...	269133	1
Bombastus	Non-ionic hydrous liniment S...	269843	1
Gehe	Non-ionic hydrous liniment S...	269845	1
Bombastus	Non-ionic hydrous liniment S...	269847	2
Noweda	Non-ionic hydrous liniment S...	269847	1
Noweda	Non-ionic hydrous liniment S...	269848	1
Noweda	Non-ionic hydrous liniment S...	271243	1
Bombastus	Non-ionic hydrous liniment S...	271245	2
Phönix Gotha	Non-ionic hydrous liniment S...	270007	2
Bombastus	Non-ionic hydrous liniment S...	271249	1
Phoenix, 12.09.2012	Non-ionic hydrous liniment S...	271302	1
Phönix Gotha	Non-ionic hydrous liniment S...	273375	2
Gehe	Non-ionic hydrous liniment S...	271305	1
Bombastus	Non-ionic hydrous liniment S...	272614	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Non-ionic hydrous liniment S...	272611	2
Bombastus	Non-ionic hydrous liniment S...	273375	1
Gehe	Non-ionic hydrous liniment S...	274014	1
Bombastus	Non-ionic hydrous liniment S...	274016	2
Apomix/Phönix	Non-ionic hydrous liniment S...	9A290513	1
Apomix 001	Non-ionic hydrous liniment S...	10A150512	2
Apomix/ Phönix	Non-ionic hydrous liniment S...	11A010713	1
apomix/Phönix	Non-ionic hydrous liniment S...	17A030912	2
Caelo	Non-ionic hydrous liniment S...	16113901	6
Apomix	Non-ionic hydrous liniment S...	20A051113	1
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	17A230913	1
Apomix / Phönix	Non-ionic hydrous liniment S...	20A051113	1
Apomix	Non-ionic hydrous liniment S...	24A101212	1
Pharmachem / Anzag	Non-ionic hydrous liniment S...	250512	1
Bombastus	Non-ionic hydrous liniment S...	270673	2
Bombastus	Non-ionic hydrous liniment S...	275626	1
Sanacorp	Non-ionic hydrous liniment S...	274820	1
Bombastus	Non-ionic hydrous liniment S...	275785	3
Bombastus	Non-ionic hydrous liniment S...	275788	6
Bombastus	Non-ionic hydrous liniment S...	276525	6
Phönix 05.04.2013	Non-ionic hydrous liniment S...	275787	1
Bombastus	Non-ionic hydrous liniment S...	278256	1
Bombastus	Non-ionic hydrous liniment S...	278260	5
Bombastus	Non-ionic hydrous liniment S...	279496	3
Bombastus	Non-ionic hydrous liniment S...	279502	6
Phönix Gotha	Non-ionic hydrous liniment S...	281559	2
Bombastus	Non-ionic hydrous liniment S...	280876	2
Caelo	Non-ionic hydrous liniment S...	281564	2
Bom/Kehr	Non-ionic hydrous liniment S...	282332	1
Bombastus	Non-ionic hydrous liniment S...	282332	4
Bombastus	Hydrated liniment SR/DAC (NR...	266553	1
Gehe	Hydrated liniment SR/DAC (NR...	269845	1
Bombastus	Hydrated liniment SR/DAC (NR...	268572	1
Bombastus	Hydrated liniment SR/DAC (NR...	268577	2
Bombastus	Hydrated liniment SR/DAC (NR...	269848	1
Bombastus	Hydrated liniment SR/DAC (NR...	270479	1
Bombastus	Hydrated liniment SR/DAC (NR...	270480	1
Bombastus	Hydrated liniment SR/DAC (NR...	272646	1
Bombastus	Non-ionic hydrous liniment S...	272615	4
Bombastus	Hydrated liniment SR/DAC (NR...	273772	2
Phönix Gotha	Hydrated liniment SR/DAC (NR...	274014	1
Bombastus	Hydrated liniment SR/DAC (NR...	273774	3
Bombastus	Non-ionic hydrous liniment S...	274015	1
Noweda	Non-ionic hydrous liniment S...	274016	1
Bombastus	Non-ionic hydrous liniment S...	274014	2
Pharmachem Gmbh	Hydrated liniment SR/DAC (NR...	20614	2
Angag Leyh	Hydrated liniment SR/DAC (NR...	4580715	1
apomix/Phönix	Hydrated liniment SR/DAC (NR...	5A210113	1
Pharmachem/Phönix	Hydrated liniment SR/DAC (NR...	21013	1
PKH Halle	Non-ionic hydrous liniment S...	6A120312	1
Apomix	Hydrated liniment SR/DAC (NR...	10A260213	1
Anzag	Non-ionic hydrous liniment S...	91012	1
Sanacorp/ Leyh Pharma	Hydrated liniment SR/DAC (NR...	115051	1
Apomix	Hydrated liniment SR/DAC (NR...	6A2000114	1
Leyh Pharma GmbH	Hydrated liniment SR/DAC (NR...	110531	1
Leyh / Anzag	Hydrated liniment SR/DAC (NR...	122531	4
Leyh Pharma	Hydrated liniment SR/DAC (NR...	122141	1
Sanacorp/ Leyh Pharma	Hydrated liniment SR/DAC (NR...	122141	1
leyh Pharma/ Sanacorp	Hydrated liniment SR/DAC (NR...	124941	1
Leyh Pharma	Hydrated liniment SR/DAC (NR...	123621	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Leyh/Gehe	Hydrated liniment SR/DAC (NR...	131121	1
Leyh	Hydrated liniment SR/DAC (NR...	124941	1
apomix	Hydrated liniment SR/DAC (NR...	12A150212	1
Sanacorp/ Leyh Pharma	Hydrated liniment SR/DAC (NR...	123011	1
Leyh Pharma	Hydrated liniment SR/DAC (NR...	124941	1
pkh, Phönix	Hydrated liniment SR/DAC (NR...	14B170214	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	15A050313	1
Apomix/Phönix	Hydrated liniment SR/DAC (NR...	15A260813	2
Apomix	Hydrated liniment SR/DAC (NR...	16A070312	2
Apomix	Non-ionic hydrous liniment S...	17A030912	1
apomix	Hydrated liniment SR/DAC (NR...	23A020413	1
PKH GmbH/ Phönix	Hydrated liniment SR/DAC (NR...	23A020413	1
Pharmachem/ANZAG/27.03...	Hydrated liniment SR/DAC (NR...	210113	1
Apomix	Hydrated liniment SR/DAC (NR...	20B180313	1
Apomix	Hydrated liniment SR/DAC (NR...	17B050313	1
Bombastus	Hydrated liniment SR/DAC (NR...	274772	3
Bombastus	Non-ionic hydrous liniment S...	274869	1
Noweda	Non-ionic hydrous liniment S...	275786	1
Bombastus	Non-ionic hydrous liniment S...	275787	2
Phönix Gotha	Hydrated liniment SR/DAC (NR...	275789	1
Bombastus	Hydrated liniment SR/DAC (NR...	275792	1
Phönix Gotha	Non-ionic hydrous liniment S...	275792	1
Bombastus	Hydrated liniment SR/DAC (NR...	275869	1
Bombastus	Hydrated liniment SR/DAC (NR...	275864	1
Bombastus	Hydrated liniment SR/DAC (NR...	277781	2
Bombastus	Hydrated liniment SR/DAC (NR...	277782	2
Bombastus	Hydrated liniment SR/DAC (NR...	279398	2
Bombastus	Hydrated liniment SR/DAC (NR...	279400	4
Phönix Gotha	Hydrated liniment SR/DAC (NR...	279496	1
Apomix	Hydrated liniment SR/DAC (NR...	27A230413	1
apomix	Hydrated liniment SR/DAC (NR...	27B230413	1
Bombastus	Hydrated liniment SR/DAC (NR...	280805	2
Bombastus	Hydrated liniment SR/DAC (NR...	281642	3
Bombastus	Hydrated liniment SR/DAC (NR...	281562	1
Bombstus/ Sanacorp	Hydrated liniment SR/DAC (NR...	281643	2
Bombastus	Hydrated liniment SR/DAC (NR...	281645	2
Bombastus	Hydrated liniment SR/DAC (NR...	281646	2
Bombastus	Hydrated liniment SR/DAC (NR...	282884	1
Bombastus	Hydrated liniment SR/DAC (NR...	284338	1
Apomix/Kehr	Hydrated liniment SR/DAC (NR...	29A060513	1
pkh, Phönix	Hydrated liniment SR/DAC (NR...	32A120515	1
Anzag	Hydrated liniment SR/DAC (NR...	30B130513	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	40A170613	1
apomix	Hydrated liniment SR/DAC (NR...	33A210513	1
Apomix / Phönix	Hydrated liniment SR/DAC (NR...	38A040613	1
Apomix	Hydrated liniment SR/DAC (NR...	36B270513	1
Euro OTC	Hydrated liniment SR/DAC (NR...	33A210513	1
Apomix	Hydrated liniment SR/DAC (NR...	42B240613	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	43A010713	1
Apomix	Hydrated liniment SR/DAC (NR...	43A010713	1
apomix	Hydrated liniment SR/DAC (NR...	46B150713	1
Caelo	Hydrated liniment SR/DAC (NR...	47B020712	1
Apomix	Hydrated liniment SR/DAC (NR...	49B170713	1
Apomix/Phoenix	Hydrated liniment SR/DAC (NR...	51A230713	1
Apomix	Hydrated liniment SR/DAC (NR...	51A230713	2
Apomix Anzag	Hydrated liniment SR/DAC (NR...	51B170712	2
PKH Halle	Hydrated liniment SR/DAC (NR...	56A120813	4
PHK GmbH	Hydrated liniment SR/DAC (NR...	66A160913	1
Apomix	Hydrated liniment SR/DAC (NR...	62B050913	1
apomix	Hydrated liniment SR/DAC (NR...	56A240712	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix/Phönix	Hydrated liniment SR/DAC (NR...	68A300913	1
Apomix	Hydrated liniment SR/DAC (NR...	70A071013	2
PHK GmbH	Hydrated liniment SR/DAC (NR...	71A291013	2
apomix	Hydrated liniment SR/DAC (NR...	72B041113	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	71A291013	2
Caelo	Hydrated liniment SR/DAC (NR...	74B081012	2
Phoenix/Apomix	Hydrated liniment SR/DAC (NR...	71A291013	1
Gehe/Apomix	Hydrated liniment SR/DAC (NR...	77A251113	1
Apomix	Hydrated liniment SR/DAC (NR...	76B1611012	1
apomix	Hydrated liniment SR/DAC (NR...	79A091213	1
Apomix	Hydrated liniment SR/DAC (NR...	79A091213	2
Apomix	Hydrated liniment SR/DAC (NR...	83A071112	1
Apomix/Gehe	Hydrated liniment SR/DAC (NR...	85A281112	1
apomix	Hydrated liniment SR/DAC (NR...	88A061211	1
Apomix	Hydrated liniment SR/DAC (NR...	90A171212	1
PKH GmbH / Phönix	Non-ionic hydrous liniment S...	11B270515	1
PharmaChem	Non-ionic hydrous liniment S...	140814	1
Caelo	Non-ionic hydrous liniment S...	15014601	8
Caelo	Non-ionic hydrous liniment S...	15096101	11
Caelo	Non-ionic hydrous liniment S...	15192301	7
Caelo	Non-ionic hydrous liniment S...	16021601	13
Bombastus	Non-ionic hydrous liniment S...	285025	2
Bombastus	Hydrated liniment SR/DAC (NR...	285740	1
Bombastus	Non-ionic hydrous liniment S...	288451	5
Bombastus	Hydrated liniment SR/DAC (NR...	288536	2
Bombastus	Hydrated liniment SR/DAC (NR...	289863	3
Bombastus	Non-ionic hydrous liniment S...	290612	3
Bombastus	Non-ionic hydrous liniment S...	291286	4
Bombastus	Hydrated liniment SR/DAC (NR...	292493	4
Bombastus	Hydrated liniment SR/DAC (NR...	294435	2
Bombastus	Non-ionic hydrous liniment S...	295069	3
Bombastus	Non-ionic hydrous liniment S...	295073	1
Bombastus	Non-ionic hydrous liniment S...	267329	1
Bombastus	Non-ionic hydrous liniment S...	269131	1
Anzag	Non-ionic hydrous liniment S...	274016	1
Bombastus	Non-ionic hydrous liniment S...	271236	1
Bombastus	Non-ionic hydrous liniment S...	273371	1
Bombastus	Non-ionic hydrous liniment S...	274813	1
Bombastus	Non-ionic hydrous liniment S...	274819	2
Gehe	Non-ionic hydrous liniment S...	275792	1
Gehe	Non-ionic hydrous liniment S...	279498	2
Bombastus	Non-ionic hydrous liniment S...	287320	1
Gehe	Non-ionic hydrous liniment S...	289205	1
	Non-ionic hydrous liniment S...	288452	1
Bombastus	Hydrated liniment SR/DAC (NR...	265946	1
Bombastus	Hydrated liniment SR/DAC (NR...	273773	1
Bombastus	Hydrated liniment SR/DAC (NR...	270673	1
Sanacorp	Hydrated liniment SR/DAC (NR...	275869	1
Bombastus	Hydrated liniment SR/DAC (NR...	279399	1
Apomix	Hydrated liniment SR/DAC (NR...	77A251113	1
PKH GmbH/ Sanacorp	Hydrated liniment SR/DAC (NR...	78B031213	1
Caelo	Non-ionic hydrous liniment S...	275785	1
Bombastus	Non-ionic hydrous liniment S...	276523	1
Anzag	Non-ionic hydrous liniment S...	275786	1
Bombastus	Non-ionic hydrous liniment S...	278255	1
Sanacorp	Non-ionic hydrous liniment S...	279500	1
Bombastus/Sanac	Non-ionic hydrous liniment S...	279498	1
Bombastus	Non-ionic hydrous liniment S...	279500	3
Bombastus	Non-ionic hydrous liniment S...	279501	4
Bombastus	Non-ionic hydrous liniment S...	281563	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Non-ionic hydrous liniment S...	282946	3
Bombastus	Non-ionic hydrous liniment S...	282334	4
Bombastus	Non-ionic hydrous liniment S...	304936	2
Gehe Kassel	Non-ionic hydrous liniment S...	15014601	2
Apomix / Gehe	Non-ionic hydrous liniment S...	1A100117	1
Apomix	Non-ionic hydrous liniment S...	1A100117	1
apomix / Anzag	Non-ionic hydrous liniment S...	1A200115	1
Apomix	Non-ionic hydrous liniment S...	2A020216	4
Apomix	Non-ionic hydrous liniment S...	2A100215	2
GEHE	Non-ionic hydrous liniment S...	2A170118	1
apomix/Noweda	Non-ionic hydrous liniment S...	2A170118	1
Apomix/Noweda	Non-ionic hydrous liniment S...	3B180214	1
Apomix / Phönix	Non-ionic hydrous liniment S...	4A010316	1
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	4A250214	1
Apomix	Non-ionic hydrous liniment S...	4A250214	2
Caelo	Non-ionic hydrous liniment S...	5A110315	1
Apomix	Non-ionic hydrous liniment S...	5A110315	1
Apomix	Non-ionic hydrous liniment S...	5A120318	1
apomix/Phönix	Non-ionic hydrous liniment S...	5B170414	2
Caelo	Non-ionic hydrous liniment S...	6A290317	1
apomix	Non-ionic hydrous liniment S...	6A310316	1
Apomix	Non-ionic hydrous liniment S...	6A310316	2
Apomix/Noweda	Non-ionic hydrous liniment S...	6B180315	1
Noweda/Apomix	Non-ionic hydrous liniment S...	6B180315	1
apomix	Non-ionic hydrous liniment S...	7A170414	1
Apomix	Non-ionic hydrous liniment S...	7A170414	1
PKH GmbH	Non-ionic hydrous liniment S...	7A300316	1
apomix/Phönix	Non-ionic hydrous liniment S...	8A080415	1
Noweda	Non-ionic hydrous liniment S...	8A080415	1
apomix	Non-ionic hydrous liniment S...	8A100516	1
apomix	Non-ionic hydrous liniment S...	9A190516	1
apomix PKH GmbH Halle	Non-ionic hydrous liniment S...	9A190516	1
PKH/Phönix	Non-ionic hydrous liniment S...	9A190516	1
Apomix	Non-ionic hydrous liniment S...	9A190516	1
Caelo	Non-ionic hydrous liniment S...	102600630	1
AHD 09.01.15/apomix	Non-ionic hydrous liniment S...	A150914	1
apomix	Non-ionic hydrous liniment S...	10A110515	1
Caelo	Non-ionic hydrous liniment S...	10A110515	1
Caelo	Non-ionic hydrous liniment S...	10260682	1
Apomix / Phönix	Non-ionic hydrous liniment S...	10A150514	1
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	10A150514	1
Caelo	Non-ionic hydrous liniment S...	10A310516	1
Apomix/ Anzag	Non-ionic hydrous liniment S...	10A310516	1
Noweda	Non-ionic hydrous liniment S...	10A150514	1
Apomix	Non-ionic hydrous liniment S...	11B270515	1
apomix	Non-ionic hydrous liniment S...	12A100615	1
Noweda/apomix	Non-ionic hydrous liniment S...	12A100615	1
PKH Halle/Phönix	Non-ionic hydrous liniment S...	12A200616	1
Apomix	Non-ionic hydrous liniment S...	12A200616	1
Apomix	Non-ionic hydrous liniment S...	12A270514	1
Apomix/Noweda	Non-ionic hydrous liniment S...	13A040614	1
apomix	Non-ionic hydrous liniment S...	13A040614	1
Caelo	Non-ionic hydrous liniment S...	13A260617	1
Caelo	Non-ionic hydrous liniment S...	14024102	1
Caelo	Non-ionic hydrous liniment S...	14092015A	1
Caelo	Non-ionic hydrous liniment S...	14134601	13
Caelo	Non-ionic hydrous liniment S...	14173601	3
Caesar/Noweda	Non-ionic hydrous liniment S...	14173601	1
Caelo	Non-ionic hydrous liniment S...	14209501	9
Caelo	Non-ionic hydrous liniment S...	14276001	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Non-ionic hydrous liniment S...	14328402	3
Caelo	Non-ionic hydrous liniment S...	14350703	2
Fagron	Non-ionic hydrous liniment S...	14350703	1
Alliance health	Non-ionic hydrous liniment S...	14350703	2
	Non-ionic hydrous liniment S...	14350703	1
GEHE	Non-ionic hydrous liniment S...	14A260617	1
apomix	Non-ionic hydrous liniment S...	14A300615	1
Caelo	Non-ionic hydrous liniment S...	15002302	5
Sanacorp	Non-ionic hydrous liniment S...	15002302	1
Sanacorp	Non-ionic hydrous liniment S...	15014601	1
Caelo	Non-ionic hydrous liniment S...	15096001	1
Cealo/Noweda	Non-ionic hydrous liniment S...	15096101	1
Noweda 27.08.15	Non-ionic hydrous liniment S...	15192301	1
Sanacorp 21.10.15 und...	Non-ionic hydrous liniment S...	15192301	1
Euro OTC	Non-ionic hydrous liniment S...	15192401	1
Bombastus	Non-ionic hydrous liniment S...	15192502	1
Caelo	Non-ionic hydrous liniment S...	15279702	1
Cael/Phönix	Non-ionic hydrous liniment S...	15326702	1
Caesar/Noweda	Non-ionic hydrous liniment S...	15326703	1
Leyh Pharma	Non-ionic hydrous liniment S...	153811	1
Caelo	Non-ionic hydrous liniment S...	15420901	2
Noweda	Non-ionic hydrous liniment S...	15420902	1
Caelo	Non-ionic hydrous liniment S...	15420902	3
Noweda	Non-ionic hydrous liniment S...	15421002	1
Caelo	Non-ionic hydrous liniment S...	15421002	3
wepa/Sanaco	Non-ionic hydrous liniment S...	1542102	1
Eigene Herstellung	Non-ionic hydrous liniment S...	1546	1
Apomix / Phönix	Non-ionic hydrous liniment S...	15A140714	1
Eigene Herstellung	Non-ionic hydrous liniment S...	1567	1
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	15A140714	1
apomix / Anzag	Non-ionic hydrous liniment S...	15A140714	1
apomix	Non-ionic hydrous liniment S...	15A150715	2
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	15A150715	1
Caelo	Non-ionic hydrous liniment S...	16113801	8
Celo	Non-ionic hydrous liniment S...	16113801	1
Caelo	Non-ionic hydrous liniment S...	16113902	1
Caelo	Non-ionic hydrous liniment S...	16258401	5
Caelo	Non-ionic hydrous liniment S...	16258501	6
hochhaus-Apotheke	Non-ionic hydrous liniment S...	16318801	1
Caelo	Non-ionic hydrous liniment S...	16318802	4
Noweda	Non-ionic hydrous liniment S...	16318802	1
Apomix	Non-ionic hydrous liniment S...	16A040814	1
apomix	Non-ionic hydrous liniment S...	16A050815	2
GEHE	Non-ionic hydrous liniment S...	16A020817	1
Apomix	Non-ionic hydrous liniment S...	16A050815	2
Caelo	Non-ionic hydrous liniment S...	16A050815	2
Apomix / Phönix	Non-ionic hydrous liniment S...	16A050815	1
Bombastus	Non-ionic hydrous liniment S...	17000577	1
Bombastus	Non-ionic hydrous liniment S...	17001403	1
Bombastus	Non-ionic hydrous liniment S...	17001411	1
Bombastus	Non-ionic hydrous liniment S...	17001412	1
Cealo/Noweda	Non-ionic hydrous liniment S...	17027801	1
phönix	Non-ionic hydrous liniment S...	17027801	1
Fagron	Non-ionic hydrous liniment S...	17027901	2
Caelo	Non-ionic hydrous liniment S...	17064401	3
Gehe	Non-ionic hydrous liniment S...	17064404	1
Caelo	Non-ionic hydrous liniment S...	17144401	1
Caelo	Non-ionic hydrous liniment S...	17144502	1
Caelo	Non-ionic hydrous liniment S...	17208901	1
Caelo	Non-ionic hydrous liniment S...	17312501	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
GEHE	Non-ionic hydrous liniment S...	17312501	1
Caelo	Non-ionic hydrous liniment S...	17312601	2
Caelo	Non-ionic hydrous liniment S...	17313201	1
Caelo	Non-ionic hydrous liniment S...	17A140817	1
apomix	Non-ionic hydrous liniment S...	17A220816	1
Bombastus	Non-ionic hydrous liniment S...	18000057	2
Bombastus	Non-ionic hydrous liniment S...	18000058	2
Bombastus	Non-ionic hydrous liniment S...	18000059	2
Bombastus	Non-ionic hydrous liniment S...	18001767	2
Bombastus	Non-ionic hydrous liniment S...	18002332	1
GEHE	Non-ionic hydrous liniment S...	18050701	1
apomix / Gehe	Non-ionic hydrous liniment S...	18B010916	1
Apomix	Non-ionic hydrous liniment S...	18A200814	1
apomix/Noweda	Non-ionic hydrous liniment S...	19A010916	1
Apomix	Non-ionic hydrous liniment S...	19A010916	1
AHD 09.01.15/apomix	Non-ionic hydrous liniment S...	20A150914	1
apomix	Non-ionic hydrous liniment S...	20A280915	1
Apomix	Non-ionic hydrous liniment S...	20A280915	1
GEHE	Non-ionic hydrous liniment S...	21A250917	1
PKH GmbH	Non-ionic hydrous liniment S...	21A241016	1
apomix / Anzag	Non-ionic hydrous liniment S...	22A101114	1
Apomix	Non-ionic hydrous liniment S...	22A101114	2
apomix/Phönix	Non-ionic hydrous liniment S...	22A101114	1
AHD 09.01.15/apomix	Non-ionic hydrous liniment S...	22A101114	1
Noweda	Non-ionic hydrous liniment S...	22A101114	1
PKH GmbH	Non-ionic hydrous liniment S...	22A251016	1
Apomix	Non-ionic hydrous liniment S...	23A181115	2
apomix AHD	Non-ionic hydrous liniment S...	23A181115	1
Pöskneek Pharma / Gehe	Non-ionic hydrous liniment S...	23A181115	2
GEHE	Non-ionic hydrous liniment S...	24A231116	1
Caelo	Non-ionic hydrous liniment S...	24A271114	1
Apomix	Non-ionic hydrous liniment S...	25A021214	2
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	25A021214	1
Gehe	Non-ionic hydrous liniment S...	25A071217	1
Apomix/ Phönix	Non-ionic hydrous liniment S...	25A081215	1
Apomix	Non-ionic hydrous liniment S...	25A081215	1
PKH Halle/Noweda	Non-ionic hydrous liniment S...	25A081215	1
Noweda	Non-ionic hydrous liniment S...	25A191213	1
Apomix / Phönix	Non-ionic hydrous liniment S...	25A191213	3
apomix	Non-ionic hydrous liniment S...	26A281215	1
Bombastus	Non-ionic hydrous liniment S...	26A281216	1
GEHE	Non-ionic hydrous liniment S...	27A181217	1
Apomix	Non-ionic hydrous liniment S...	27A221214	2
Apomix/Gehe	Non-ionic hydrous liniment S...	27A221214	1
gehe	Non-ionic hydrous liniment S...	281560	2
Caelo	Non-ionic hydrous liniment S...	281562	1
Caelo	Non-ionic hydrous liniment S...	282946	1
Bombastus	Non-ionic hydrous liniment S...	282335	1
Bombastus	Non-ionic hydrous liniment S...	282333	4
Bombastus	Non-ionic hydrous liniment S...	285023	1
Bombastus	Non-ionic hydrous liniment S...	285024	1
Bombastu-Werke AG/Gehe	Non-ionic hydrous liniment S...	285025	1
Caelo	Non-ionic hydrous liniment S...	285025	1
Bombastus	Non-ionic hydrous liniment S...	285790	2
Bombastus	Non-ionic hydrous liniment S...	285792	2
Caelo	Non-ionic hydrous liniment S...	285791	2
Wolf-Apotheke	Non-ionic hydrous liniment S...	285793	1
gehe	Non-ionic hydrous liniment S...	285793	1
Bombastus	Non-ionic hydrous liniment S...	285795	3
Noweda	Non-ionic hydrous liniment S...	285795	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Non-ionic hydrous liniment S...	285838	1
Gehe/Apomix	Non-ionic hydrous liniment S...	286417	1
Bombastus	Non-ionic hydrous liniment S...	286420	1
Apomix	Non-ionic hydrous liniment S...	286420	1
Gehe	Non-ionic hydrous liniment S...	286421	1
Gehe	Non-ionic hydrous liniment S...	287316	1
Bombastus	Non-ionic hydrous liniment S...	287319	2
Caelo	Non-ionic hydrous liniment S...	287323	1
Noweda	Non-ionic hydrous liniment S...	287323	1
Bombastus	Non-ionic hydrous liniment S...	288386	1
Bombastus	Non-ionic hydrous liniment S...	288387	2
Gehe 18.08.14	Non-ionic hydrous liniment S...	288387	4
Gehe	Non-ionic hydrous liniment S...	288387	1
Bombastus	Non-ionic hydrous liniment S...	288389	2
Bombastus	Non-ionic hydrous liniment S...	288390	1
Bombastu-Werke AG/Gehe	Non-ionic hydrous liniment S...	288452	1
Bombastus	Non-ionic hydrous liniment S...	288452	1
Bombastus	Non-ionic hydrous liniment S...	289205	1
Bombastus	Non-ionic hydrous liniment S...	289206	3
Bombastus	Non-ionic hydrous liniment S...	289208	2
Bombastus	Non-ionic hydrous liniment S...	289207	1
Bomastus/Phönix	Non-ionic hydrous liniment S...	290610	1
Noweda	Non-ionic hydrous liniment S...	290610	2
Bombastus	Non-ionic hydrous liniment S...	291220	2
Bombastus	Non-ionic hydrous liniment S...	291221	3
Caelo	Non-ionic hydrous liniment S...	291222	3
Caelo	Non-ionic hydrous liniment S...	291286	1
Bombastus	Non-ionic hydrous liniment S...	291861	3
Bombastus	Non-ionic hydrous liniment S...	292617	1
Noweda	Non-ionic hydrous liniment S...	292618	1
Bombastus	Non-ionic hydrous liniment S...	292621	3
Caelo	Non-ionic hydrous liniment S...	292742	1
Bombastus	Non-ionic hydrous liniment S...	293574	1
Bombastus	Non-ionic hydrous liniment S...	293577	3
Bombastus	Non-ionic hydrous liniment S...	293578	2
Bombastus	Non-ionic hydrous liniment S...	293579	5
Caelo	Non-ionic hydrous liniment S...	293579	1
Caelo	Non-ionic hydrous liniment S...	294211	2
Bombastus	Non-ionic hydrous liniment S...	294212	1
Bombastus	Non-ionic hydrous liniment S...	294213	2
Bombastus	Non-ionic hydrous liniment S...	294277	2
Bombastus	Non-ionic hydrous liniment S...	295068	1
Bombastus	Non-ionic hydrous liniment S...	295070	6
Bombastus	Non-ionic hydrous liniment S...	295072	2
Bombastus	Non-ionic hydrous liniment S...	295071	1
Caelo	Non-ionic hydrous liniment S...	295072	1
Bombastus	Non-ionic hydrous liniment S...	295910	3
Caelo	Non-ionic hydrous liniment S...	295910	1
Bombastus	Non-ionic hydrous liniment S...		2
Bombastus	Non-ionic hydrous liniment S...	296107	2
Bombastus	Non-ionic hydrous liniment S...	296108	1
Bombastus	Non-ionic hydrous liniment S...	296109	2
Noweda	Non-ionic hydrous liniment S...	296109	1
Bombastus	Non-ionic hydrous liniment S...	296111	1
Bombastus	Non-ionic hydrous liniment S...	296584	7
Bombastus	Non-ionic hydrous liniment S...	296585	7
Caelo	Non-ionic hydrous liniment S...	296585	1
Bombastus	Non-ionic hydrous liniment S...	296589	3
Bombastus	Non-ionic hydrous liniment S...	297171	1
Bombastus	Non-ionic hydrous liniment S...	297172	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Non-ionic hydrous liniment S...	297174	2
Bombastus	Non-ionic hydrous liniment S...	297173	2
Rondell Apotheke	Non-ionic hydrous liniment S...	297175	1
Bombastus	Non-ionic hydrous liniment S...	297176	2
Bombastus	Non-ionic hydrous liniment S...	297178	4
Caelo	Non-ionic hydrous liniment S...	297178	1
Bombastus	Non-ionic hydrous liniment S...	298465	2
Bombastus	Non-ionic hydrous liniment S...	298466	1
Caelo	Non-ionic hydrous liniment S...	298466	1
Bombastus	Non-ionic hydrous liniment S...	298467	1
Bombastus	Non-ionic hydrous liniment S...	299132	1
Bombastus	Non-ionic hydrous liniment S...	299133	1
Bombastus	Non-ionic hydrous liniment S...	300196	2
Bombastus	Non-ionic hydrous liniment S...	300198	1
Caelo	Non-ionic hydrous liniment S...	300199	3
Bombastus	Non-ionic hydrous liniment S...	300199	1
Bombastus	Non-ionic hydrous liniment S...	300199	2
Bombastus	Non-ionic hydrous liniment S...	300858	1
Bombastus	Non-ionic hydrous liniment S...	300859	5
Bombastus	Non-ionic hydrous liniment S...	300860	3
Bombastus	Non-ionic hydrous liniment S...	300861	2
Caelo	Non-ionic hydrous liniment S...	300864	2
Bombastus	Non-ionic hydrous liniment S...	301555	1
Bombastus	Non-ionic hydrous liniment S...	301556	2
Caelo	Non-ionic hydrous liniment S...	302283	1
Bombastus	Non-ionic hydrous liniment S...	302283	1
Phoenix	Non-ionic hydrous liniment S...	302283	1
Bombastus	Non-ionic hydrous liniment S...	302284	3
Bombastus	Non-ionic hydrous liniment S...	302285	2
Caelo	Non-ionic hydrous liniment S...	302285	1
Beiersdorf/Kehr	Non-ionic hydrous liniment S...	302849	1
Bombastus	Non-ionic hydrous liniment S...	302849	2
Bombastus	Non-ionic hydrous liniment S...	302852	6
Bombastus	Non-ionic hydrous liniment S...	302887	1
Phoenix	Non-ionic hydrous liniment S...	303895	1
Bombastus	Non-ionic hydrous liniment S...	303895	2
Bombastus	Non-ionic hydrous liniment S...	304496	1
Bombastus	Non-ionic hydrous liniment S...	303896	1
Bombastus	Non-ionic hydrous liniment S...	304935	5
Bombastus	Non-ionic hydrous liniment S...	305593	2
Bombastus	Non-ionic hydrous liniment S...	306358	3
Bombastus	Non-ionic hydrous liniment S...	306359	2
Bombastus	Non-ionic hydrous liniment S...	306360	1
Caelo	Non-ionic hydrous liniment S...	307658	1
Bombastus	Non-ionic hydrous liniment S...	307681	1
Bombastus	Non-ionic hydrous liniment S...	307683	1
Caelo	Non-ionic hydrous liniment S...	307685	1
Bombastus	Non-ionic hydrous liniment S...	307685	2
Bombastus	Non-ionic hydrous liniment S...	308623	2
Bombastus	Non-ionic hydrous liniment S...	308624	2
Caelo	Non-ionic hydrous liniment S...	308627	1
Bombastus	Non-ionic hydrous liniment S...	308627	1
Bombastus	Non-ionic hydrous liniment S...	30898	1
Bombastus	Non-ionic hydrous liniment S...	309790	1
Bombastus	Non-ionic hydrous liniment S...	309795	3
Bombastus	Non-ionic hydrous liniment S...	309798	1
Bombastus	Non-ionic hydrous liniment S...	309804	4
Bombastus	Non-ionic hydrous liniment S...	310355	1
Fagron	Non-ionic hydrous liniment S...	311041	1
Bombastus	Non-ionic hydrous liniment S...	311041	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Non-ionic hydrous liniment S...	311609	2
Caelo	Non-ionic hydrous liniment S...	4063101	1
Bombastus	Non-ionic hydrous liniment S...	60091111	1
Bombastus	Non-ionic hydrous liniment S...	98434	1
Bombastus	Non-ionic hydrous liniment S...	60091112	1
Bombastus	Non-ionic hydrous liniment S...	229263	1
Bombastus	Non-ionic hydrous liniment S...	267339	1
Bombastus	Non-ionic hydrous liniment S...	268124	1
Gehe	Non-ionic hydrous liniment S...	269846	1
Bombastus	Non-ionic hydrous liniment S...	269846	1
Gehe	Non-ionic hydrous liniment S...	271303	1
Bombastus	Non-ionic hydrous liniment S...	271304	2
Bombastus	Non-ionic hydrous liniment S...	273374	1
Bombastus	Non-ionic hydrous liniment S...	273373	1
Fagron	Non-ionic hydrous liniment S...	27301	1
Apomix Anzag	Non-ionic hydrous liniment S...	6A080413	1
PKH Halle / Phoenix	Non-ionic hydrous liniment S...	4A110313	1
apomix	Non-ionic hydrous liniment S...	8A170412	1
Bombastus	Non-ionic hydrous liniment S...	7021402	1
apomix	Non-ionic hydrous liniment S...	10A1206013	1
apomix	Non-ionic hydrous liniment S...	11A010713	1
Apomix	Non-ionic hydrous liniment S...	11A010713	1
Fagron	Non-ionic hydrous liniment S...	122304	1
Caelo	Non-ionic hydrous liniment S...	122304	1
Caelo	Non-ionic hydrous liniment S...	123304	1
Fagron	Non-ionic hydrous liniment S...	123304	2
Fagron	Non-ionic hydrous liniment S...	125304	1
Apomix/Noweda	Non-ionic hydrous liniment S...	12A100713	1
Apomix Anzag	Non-ionic hydrous liniment S...	12A270514	1
Bombastus	Non-ionic hydrous liniment S...	131648	1
Caelo	Non-ionic hydrous liniment S...	14024104	1
Pharmachem / Anzag	Non-ionic hydrous liniment S...	150213	1
Phönix/Apomix	Non-ionic hydrous liniment S...	15A260813	1
Caelo	Non-ionic hydrous liniment S...	15192391	1
Apomix	Non-ionic hydrous liniment S...	15A260813	1
Apomix	Non-ionic hydrous liniment S...	15A300712	1
Phönix (PKH Halle)	Non-ionic hydrous liniment S...	16A060812	1
PKH Halle,Phönix	Non-ionic hydrous liniment S...	17A030912	1
Apomix	Non-ionic hydrous liniment S...	17A230913	1
Apomix / Phönix	Non-ionic hydrous liniment S...	17A230913	1
Caelo	Non-ionic hydrous liniment S...	2/AD1	1
apomix	Non-ionic hydrous liniment S...	20A051113	1
apomix/Phönix	Non-ionic hydrous liniment S...	20A121211	1
Selbstersteller	Non-ionic hydrous liniment S...	20140116	1
PKH Halle,Anzag	Non-ionic hydrous liniment S...	20A221012	1
Apomix / Phönix	Non-ionic hydrous liniment S...	23A111213	1
apomix	Non-ionic hydrous liniment S...	24A101212	1
Caelo	Non-ionic hydrous liniment S...	21A201113	1
apomix	Non-ionic hydrous liniment S...	20A221012	1
Fagron	Non-ionic hydrous liniment S...	257210	1
Bombastus	Non-ionic hydrous liniment S...	274815	1
Fagron	Non-ionic hydrous liniment S...	258210	1
Noweda	Non-ionic hydrous liniment S...	274820	2
Bombastus	Non-ionic hydrous liniment S...	274820	1
Phoenix am 07.03.2013	Non-ionic hydrous liniment S...	275785	2
Bombastus	Non-ionic hydrous liniment S...	275786	4
Caelo	Non-ionic hydrous liniment S...	275787	1
Anzag	Non-ionic hydrous liniment S...	275788	1
Gehe	Non-ionic hydrous liniment S...	275791	1
Bombastus	Non-ionic hydrous liniment S...	275789	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Non-ionic hydrous liniment S...	275792	1
Gehe	Non-ionic hydrous liniment S...	276524	1
Bombastus	Non-ionic hydrous liniment S...	276524	1
Bombastus	Non-ionic hydrous liniment S...	276526	2
Gehe	Non-ionic hydrous liniment S...	278253	2
Phönix Gotha	Non-ionic hydrous liniment S...	278255	1
Bombastus	Non-ionic hydrous liniment S...	278258	4
Bombastus	Non-ionic hydrous liniment S...	278259	6
Caesar & Loretz/ Noweda	Non-ionic hydrous liniment S...	278260	1
Caelo	Non-ionic hydrous liniment S...	278260	1
Bombastus	Non-ionic hydrous liniment S...	279497	1
Bombastus	Non-ionic hydrous liniment S...	279498	2
Bombastus	Non-ionic hydrous liniment S...	279499	3
Sanacorp	Non-ionic hydrous liniment S...	279501	1
post apotheke	Non-ionic hydrous liniment S...	279502	2
Gehe	Non-ionic hydrous liniment S...	280873	1
Bombastus	Non-ionic hydrous liniment S...	280877	6
Gehe	Non-ionic hydrous liniment S...	280877	1
Bombastus	Non-ionic hydrous liniment S...	281558	3
Bombastus	Non-ionic hydrous liniment S...	281560	1
Bombastus	Non-ionic hydrous liniment S...	281561	2
Bombastus	Non-ionic hydrous liniment S...	281562	2
/Noweda	Non-ionic hydrous liniment S...	281562	1
Noweda	Non-ionic hydrous liniment S...	281562	1
Bombastus	Non-ionic hydrous liniment S...	281564	4
Gehe	Non-ionic hydrous liniment S...	281563	1
Bombastu-Werke AG/Gehe	Non-ionic hydrous liniment S...	281564	1
Caelo	Non-ionic hydrous liniment S...	282333	1
Gehe	Non-ionic hydrous liniment S...	282335	1
Bombastu-Werke AG/Gehe	Non-ionic hydrous liniment S...	282947	1
gehe	Non-ionic hydrous liniment S...	282947	1
Phönix Gotha	Non-ionic hydrous liniment S...	282946	1
Bombastus	Non-ionic hydrous liniment S...	284233	1
Bombastus	Non-ionic hydrous liniment S...	284231	1
Bombastus	Non-ionic hydrous liniment S...	284234	1
Bombastus	Non-ionic hydrous liniment S...	284236	1
Gehe	Non-ionic hydrous liniment S...	285024	1
Bombastus	Non-ionic hydrous liniment S...	284237	1
Caelo	Non-ionic hydrous liniment S...	300863	1
Bombastus	Non-ionic hydrous liniment S...	600991111	1
Pharmachem/Noweda	Hydrated liniment SR/DAC (NR...	10916	1
apomix	Hydrated liniment SR/DAC (NR...	2A120115	1
Caelo	Hydrated liniment SR/DAC (NR...	20614	1
Pharmachem GmbH/Phönix	Hydrated liniment SR/DAC (NR...	30815	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	3A070114	1
Apomix /Alliance	Hydrated liniment SR/DAC (NR...	3A070114	1
apomix	Hydrated liniment SR/DAC (NR...	3B130115	1
Bombastus	Hydrated liniment SR/DAC (NR...	4A070114	1
Apomix	Hydrated liniment SR/DAC (NR...	6A200114	1
Apomix/ Galenus Apotheke	Hydrated liniment SR/DAC (NR...	7A250116	1
Apomix	Hydrated liniment SR/DAC (NR...	9A270114	1
apomix/Phönix	Hydrated liniment SR/DAC (NR...	10A080216	1
Euro Pharmachem/Phönix	Hydrated liniment SR/DAC (NR...	110315	1
apomix	Hydrated liniment SR/DAC (NR...	11A070217	1
apomix/Anzag	Hydrated liniment SR/DAC (NR...	11a270217	1
Apomix	Hydrated liniment SR/DAC (NR...	11A270217	1
Apomix / Phönix	Hydrated liniment SR/DAC (NR...	12B110214	1
GEHE/LEYH Parma/09.04.2...	Hydrated liniment SR/DAC (NR...	134211	2
apomix	Hydrated liniment SR/DAC (NR...	13A110214	2
PKH / Phoenix Leipzig	Hydrated liniment SR/DAC (NR...	14A040315	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix	Hydrated liniment SR/DAC (NR...	14B220216	1
Leyh Pharma	Hydrated liniment SR/DAC (NR...	151621	1
Pharmachem	Hydrated liniment SR/DAC (NR...	160115	1
Fagron	Hydrated liniment SR/DAC (NR...	16182403	1
Leyhs	Hydrated liniment SR/DAC (NR...	162231	1
Leyh Pharma/Noweda	Hydrated liniment SR/DAC (NR...	164021	1
Leyhs	Hydrated liniment SR/DAC (NR...	163541	1
apomix/Noweda	Hydrated liniment SR/DAC (NR...	16A200317	1
Bombastus	Hydrated liniment SR/DAC (NR...	17000558	1
Caelo	Hydrated liniment SR/DAC (NR...	17022202	1
Bombastus	Hydrated liniment SR/DAC (NR...	171229WL/310644	1
Bombastus	Hydrated liniment SR/DAC (NR...	171207WL/310644	1
Bombastus	Hydrated liniment SR/DAC (NR...	18001657	1
Bombastus	Hydrated liniment SR/DAC (NR...	18001660	1
Bombastus	Hydrated liniment SR/DAC (NR...	180305WL/307264	1
apomix	Hydrated liniment SR/DAC (NR...	18B030417	1
apomix	Hydrated liniment SR/DAC (NR...	19A070316	1
apomix PKH GmbH Halle	Hydrated liniment SR/DAC (NR...	19B100314	1
Apomix	Hydrated liniment SR/DAC (NR...	19B040417	1
Gehe/Apomix	Hydrated liniment SR/DAC (NR...	20A170315	1
apomix/Phönix	Hydrated liniment SR/DAC (NR...	24A240417	1
PKH Halle	Hydrated liniment SR/DAC (NR...	24A310314	1
sanacorp	Hydrated liniment SR/DAC (NR...	24A310314	1
AMH Niemann/Noweda	Hydrated liniment SR/DAC (NR...	26A130415	1
apomix/Gehe	Hydrated liniment SR/DAC (NR...	27A050416	6
Bombastus	Hydrated liniment SR/DAC (NR...	287323	1
Caelo	Hydrated liniment SR/DAC (NR...	287410	1
Bombastus	Hydrated liniment SR/DAC (NR...	288537	1
Fagron	Hydrated liniment SR/DAC (NR...	289863	1
Bombastus	Hydrated liniment SR/DAC (NR...	289864	2
Bombastus	Hydrated liniment SR/DAC (NR...	289862	2
apomix PKH GmbH Halle	Hydrated liniment SR/DAC (NR...	28A060416	1
Bombastus	Hydrated liniment SR/DAC (NR...	291867	2
Bombastus	Hydrated liniment SR/DAC (NR...	293591	1
Fagron	Hydrated liniment SR/DAC (NR...	294434	1
Bombastus	Hydrated liniment SR/DAC (NR...	295853	1
Bombastus	Hydrated liniment SR/DAC (NR...	295852	1
Bombastus	Hydrated liniment SR/DAC (NR...	297053	1
Bombastus	Hydrated liniment SR/DAC (NR...	297054	1
Bombastus	Hydrated liniment SR/DAC (NR...	298572	1
Euro OTC	Hydrated liniment SR/DAC (NR...	300968	1
Bombastus	Hydrated liniment SR/DAC (NR...	301992	2
Bombastus	Hydrated liniment SR/DAC (NR...	301995	1
Bombastus	Hydrated liniment SR/DAC (NR...	301561	1
Bombastus	Hydrated liniment SR/DAC (NR...	301993	1
Bombastus	Hydrated liniment SR/DAC (NR...	306362	1
Noweda	Hydrated liniment SR/DAC (NR...	307263	1
Bombastus	Hydrated liniment SR/DAC (NR...	307680	1
Bombastus	Hydrated liniment SR/DAC (NR...	308912	3
Gehe	Hydrated liniment SR/DAC (NR...	310642	1
Bombastus	Hydrated liniment SR/DAC (NR...	310643	1
Pharmachem/ Noweda	Hydrated liniment SR/DAC (NR...	310715	1
Apomix/Gehe	Hydrated liniment SR/DAC (NR...	31A040515	2
apomix/Gehe	Hydrated liniment SR/DAC (NR...	33A070514	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	33A190515	1
apomix/Holdermann	Hydrated liniment SR/DAC (NR...	32A120515	1
Anzag	Hydrated liniment SR/DAC (NR...	33A190515	1
Caelo	Hydrated liniment SR/DAC (NR...	34B170516	1
PKH GmbH	Hydrated liniment SR/DAC (NR...	35A120617	1
	Hydrated liniment SR/DAC (NR...	37B190514	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix /Alliance	Hydrated liniment SR/DAC (NR...	41A260514	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	42A020614	1
Bombastus	Hydrated liniment SR/DAC (NR...	4114E-03266	1
apomix	Hydrated liniment SR/DAC (NR...	42A020614	1
Euro OTC	Hydrated liniment SR/DAC (NR...	42B160615	1
PKH Halle / Phoenix	Hydrated liniment SR/DAC (NR...	44A220615	1
Apomix / Phönix	Hydrated liniment SR/DAC (NR...	47B170614	1
Apomix	Hydrated liniment SR/DAC (NR...	4576872	1
Caelo	Hydrated liniment SR/DAC (NR...	52B010817	1
apomix/Gehe	Hydrated liniment SR/DAC (NR...	46A280616	2
Raben-Apotheke Mildenau	Hydrated liniment SR/DAC (NR...	49A230614	1
Apomix/Anzag	Hydrated liniment SR/DAC (NR...	56B100815	1
apomix/Gehe	Hydrated liniment SR/DAC (NR...	60A150816	6
apomix	Hydrated liniment SR/DAC (NR...	61A010915	1
apomix	Hydrated liniment SR/DAC (NR...	60A150816	1
apomix	Hydrated liniment SR/DAC (NR...	62A110814	1
Apomix	Hydrated liniment SR/DAC (NR...	61A110814	1
Raben-Apotheke Mildenau	Hydrated liniment SR/DAC (NR...	63A120814	2
aopmix/Phoenix	Hydrated liniment SR/DAC (NR...	62A110814	1
Apomix	Hydrated liniment SR/DAC (NR...	63B050916	1
Apomix	Hydrated liniment SR/DAC (NR...	64A050916	1
apomix / Anzag	Hydrated liniment SR/DAC (NR...	64A050916	1
apomix/Holdermann	Hydrated liniment SR/DAC (NR...	65A060916	1
Apomix/Sana	Hydrated liniment SR/DAC (NR...	66A210915	1
apomix	Hydrated liniment SR/DAC (NR...	66A210915	1
Apomix	Hydrated liniment SR/DAC (NR...	67A010914	1
apomix	Hydrated liniment SR/DAC (NR...	67A010914	1
Euro OTC	Hydrated liniment SR/DAC (NR...	71A051015	2
PKH Halle	Hydrated liniment SR/DAC (NR...	72A101016	1
apomix/Noweda	Hydrated liniment SR/DAC (NR...	68A301017	1
Apomix	Hydrated liniment SR/DAC (NR...	72A220914	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	75A131113	1
Caelo	Hydrated liniment SR/DAC (NR...	67A12092916	1
apomix	Hydrated liniment SR/DAC (NR...	73B131014	1
apomix/Holdermann	Hydrated liniment SR/DAC (NR...	76A031115	1
Apomix PKH GmbH	Hydrated liniment SR/DAC (NR...	79011215	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	78A161115	1
apomix/Noweda	Hydrated liniment SR/DAC (NR...	81A211116	6
apomix/Phönix	Hydrated liniment SR/DAC (NR...	81A211116	1
apomix/Gehe	Hydrated liniment SR/DAC (NR...	81A251114	1
AMH Niemann/Noweda	Hydrated liniment SR/DAC (NR...	83A021214	1
Apomix/Noweda	Hydrated liniment SR/DAC (NR...	85A091214	1
PKH Halle	Hydrated liniment SR/DAC (NR...	85A281215	1
Bombastus	Hydrated liniment SR/DAC (NR...	88A161214	1
Apomix/Gehe	Hydrated liniment SR/DAC (NR...	86A161214	1
Bombastus	Hydrated liniment SR/DAC (NR...	264186	1
Phönix Gotha	Non-ionic hydrous liniment S...	270008	1
Bombastus	Hydrated liniment SR/DAC (NR...	269132	1
Phönix	Hydrated liniment SR/DAC (NR...	271895	1
Bombastus	Hydrated liniment SR/DAC (NR...	268576	1
Bombastus	Non-ionic hydrous liniment S...	272611	1
Gehe/Apomix/06.11.2013	Hydrated liniment SR/DAC (NR...	21013	1
ApoMix	Hydrated liniment SR/DAC (NR...	4A110112	1
Pharmachem/Phönix	Hydrated liniment SR/DAC (NR...	50612	1
Pharmachem/Anzag/10.07....	Hydrated liniment SR/DAC (NR...	60513	1
Apomix	Hydrated liniment SR/DAC (NR...	5A210110	1
apomix	Hydrated liniment SR/DAC (NR...	6A180112	1
Phoenix (PKH Halle)	Hydrated liniment SR/DAC (NR...	5A210113	1
Phoenix (PKH Halle)	Hydrated liniment SR/DAC (NR...	7A290113	1
GEHE/Pharmachem/20.04.2...	Hydrated liniment SR/DAC (NR...	91211	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
ApoMix	Hydrated liniment SR/DAC (NR...	12A150212	1
Apomix	Hydrated liniment SR/DAC (NR...	21A180313	2
Apomix Anzag	Hydrated liniment SR/DAC (NR...	17A030314	1
apomix 004	Hydrated liniment SR/DAC (NR...	23A270312	1
Apomix/Phönix	Hydrated liniment SR/DAC (NR...	25A220413	1
Gehe/Pharmachem /13.06....	Hydrated liniment SR/DAC (NR...	250412	1
Bombastus	Non-ionic hydrous liniment S...	269848	1
Bombastus	Hydrated liniment SR/DAC (NR...	277778	1
Bombastus	Hydrated liniment SR/DAC (NR...	278256	1
Bombastus	Hydrated liniment SR/DAC (NR...	279397	1
Caelo	Hydrated liniment SR/DAC (NR...	281646	1
Caelo	Hydrated liniment SR/DAC (NR...	279398	1
Bombastus	Hydrated liniment SR/DAC (NR...	282882	1
Apomix	Hydrated liniment SR/DAC (NR...	33A210513	2
Gehe	Hydrated liniment SR/DAC (NR...	284338	1
Apomix	Hydrated liniment SR/DAC (NR...	35A290512	1
Apomix	Hydrated liniment SR/DAC (NR...	44A120612	1
apomix	Hydrated liniment SR/DAC (NR...	48A020712	1
Apomix	Hydrated liniment SR/DAC (NR...	48A170713	1
Apomix	Hydrated liniment SR/DAC (NR...	56A120813	1
Apomix	Hydrated liniment SR/DAC (NR...	58A010812	1
apomix	Hydrated liniment SR/DAC (NR...	60A020812	1
Phönix (PKH Halle)	Hydrated liniment SR/DAC (NR...	60A020812	1
apomix	Hydrated liniment SR/DAC (NR...	65B290812	1
Caelo	Hydrated liniment SR/DAC (NR...	67b40912	1
apomix	Hydrated liniment SR/DAC (NR...	69A120912	1
Apomix	Hydrated liniment SR/DAC (NR...	69A120912	1
Apomix Phönix	Hydrated liniment SR/DAC (NR...	71A291013	1
apomix	Hydrated liniment SR/DAC (NR...	80A061112	1
Apomix	Hydrated liniment SR/DAC (NR...	72A011012	1

- 89 423 spectra from 499 *Apo-Ident* customers from a total of 26 241 batches from a further 332 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Liniment ointments* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Liniment ointments* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	640	0	34 904
Type B	0	2440	0	91 334
Type C	0	1918	11	89 309

The substance/substance group *Liniment ointments* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.0625 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 99.7541 %)
Type C	100.0000 % (> 99.6586 %)	99.4298 % (> 99.2742 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30982	31939	8.08	41.55
31169	31169	0.00	46.08
31170	31170	0.00	42.59
31438	31938	4.02	43.53
31476	32090	5.98	45.02
31827	31947	2.22	45.36
31830	31947	1.56	45.70
31938	31938	0.00	44.11
31939	31939	0.00	42.72
31943	31943	0.00	45.90
31945	31945	0.00	45.63
31947	31947	0.00	45.86
32090	32090	0.00	45.25
32305	31170	5.80	42.55
32727	31947	9.55	44.60

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Lotio Cordes®**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31075-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Lotio Cordes®

Special notes

When selecting the *Lotio Cordes®* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Lotio Cordes®*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Ichthyol	Lotio Cordes®	12D014	31075	40	1403759
Ichthyol	Lotio Cordes®	14D010	31595	20	20140826
Ichthyol	Lotio Cordes®	15D002	32487	40	20150318

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 100 spectra of 3 reference samples from the substance/substance group *Lotio Cordes*[®]. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 444 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 160 spectra from 4 reference samples of the substance/substance group *Lotio Cordes*[®].
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Ichthyol	Lotio Cordes [®]	11D029	30505	40
Ichthyol	Lotio Cordes [®]	14D010	31595 [†]	40
Ichthyol	Lotio Cordes [®]	15D021	32965	40
Ichthyol	Lotio Cordes [®]	16D011	33839	40

- 93 714 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 93 spectra from 24 *Apo-Ident* customers from 18 batches from the substance/substance group *Lotio Cordes*[®].
- These include spectra of independent samples from 15 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Ichthyol/Phönix	Lotio Cordes [®]	15D001	2
Ichthyol/Phönix	Lotio Cordes [®]	14D011	3
Ichthyol/Phönix	Lotio Cordes [®]	15D002	4
Anzag/Ichthyol-Gesellsch...	Lotio Cordes [®]	13D009	1
Ichthyol	Lotio Cordes [®]	11D030	1
Caelo	Lotio Cordes [®]	10D011	20
Alliance Health 05.06.2...	Lotio Cordes [®]	11D029	1
Caelo	Lotio Cordes [®]	12081403	1
Ichthyol	Lotio Cordes [®]	13D009	1
Ichthyol-Gesellschaft, ...	Lotio Cordes [®]	13D010	1
Ichthyol	Lotio Cordes [®]	14041402	1
Ichthyol /Noweda	Lotio Cordes [®]	14D010	2
Ichthyol/gehe	Lotio Cordes [®]	14D010	2
Ichthyol/Sanacorp	Lotio Cordes [®]	14D010	1
Ichthyol	Lotio Cordes [®]	13D010	1
Ichthyol/Phönix	Lotio Cordes [®]	14D010	2
Ichthyol/Phoenix	Lotio Cordes [®]	14D010	1
Ichthyol	Lotio Cordes [®]	14D010	2
Ichthyol	Lotio Cordes [®]	14D011	1
Caelo	Lotio Cordes [®]	14D010	1
Ichthyol/GEHE	Lotio Cordes [®]	14D011	1
Ichthyol /Noweda	Lotio Cordes [®]	15D001	1
Ichthyol/Phoenix	Lotio Cordes [®]	15D001	1
Ichthyol	Lotio Cordes [®]	15D001	2
Ichthyol-Gesellschaft	Lotio Cordes [®]	15D001	1
Ichthyol Gesellschaft/N...	Lotio Cordes [®]	15d001	1
Ichthyol/Anzag	Lotio Cordes [®]	15D002	1
Ichthyol	Lotio Cordes [®]	15D002	2
Ichthyol/Phoenix	Lotio Cordes [®]	15D002	1
Caelo	Lotio Cordes [®]	15D002	1
Ichthyol-Gesellschaft	Lotio Cordes [®]	15D002	1
Euro OTC	Lotio Cordes [®]	15D022	2
Ichthyol	Lotio Cordes [®]	16021509	1
Ichthyol	Lotio Cordes [®]	16D010	1
Ichthyol-Gesellschaft	Lotio Cordes [®]	15D022	1
Ichthyol	Lotio Cordes [®]	16D010	20
Ichthyol	Lotio Cordes [®]	27041508	1
Phoenix/Stiefel	Lotio Cordes [®]	12D014	1
Ichthyol	Lotio Cordes [®]	27041507	1
ichthyol	Lotio Cordes [®]	12D014	1
Ichthyol/gehe	Lotio Cordes [®]	13D010	1
Ichthyol/Phönix	Lotio Cordes [®]	13D009	1

- 91 259 spectra from 501 *Apo-Ident* customers from a total of 26 991 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Lotio Cordes*[®] can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Lotio Cordes*[®] and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	100	0	35 444
Type B	0	139	21	93 714
Type C	0	86	7	91 259

The substance/substance group *Lotio Cordes*[®] can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 94.0000 %)
Type B	100.0000 % (> 99.9761 %)	86.8750 % (> 85.0000 %)
Type C	100.0000 % (> 99.6587 %)	92.4731 % (> 89.2473 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31075	31075	0.00	30.90
31595	31075	3.96	28.97
32487	31075	6.91	38.77

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Macrogol 8 stereate type I**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30554-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Macrogol 8 stereate type I; Polyäthylenglycol-8-stearat

Special notes

When selecting the *Macrogol 8 stereate type I* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Macrogol 8 stereate type I*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Macrogol 8 stere...	13077501	31165	40	1407187
Caelo	Macrogol 8 stere...	13324506	31453	40	AR-15-FG-001640-01
Caelo	Macrogol 8 stere...	172186	33717	40	20170906*
Caelo	Macrogol 8 stere...	173095	33900	40	20180111*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Macrogol 8 stereate type I*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Macrogol 8 stereate type I*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Macrogol 8 stereate type I	11137510	30554	40
Caelo	Macrogol 8 stereate type I	162302	32339	40
Caelo	Macrogol 8 stereate type I	162302	32701	40

- 93 754 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 49 spectra from 25 *Apo-Ident* customers from 24 batches from the substance/substance group *Macrogol 8 stearate type I*.
- These include spectra of independent samples from 22 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macrogol 8 stearate type I	13324504	4
sana	Macrogol 8 stearate type I	13077505	1
Caelo	Macrogol 8 stearate type I	13304501	1
Caelo	Macrogol 8 stearate type I	13324501	2
Phönix	Macrogol 8 stearate type I	13324501	1
Caelo	Macrogol 8 stearate type I	13324506	3
Fagron	Macrogol 8 stearate type I	13H08-N04	1
Caelo	Macrogol 8 stearate type I	14074603	5
Fagron	Macrogol 8 stearate type I	1515Q-03130	1
Fagron	Macrogol 8 stearate type I	14A17-B08-290544	1
Caelo	Macrogol 8 stearate type I	15249601	3
Caeol	Macrogol 8 stearate type I	15249601	1
Caelo	Macrogol 8 stearate type I	15377301	4
Caelo	Macrogol 8 stearate type I	16063202	1
Gehe	Macrogol 8 stearate type I	16063202	1
Caelo	Macrogol 8 stearate type I	16063205	2
Caelo	Macrogol 8 stearate type I	17218602	2
Caelo	Macrogol 8 stearate type I	T30911	1
Noweda	Macrogol 8 stearate type I	11300607	1
Fagron	Macrogol 8 stearate type I	T32431	1
Fagron	Macrogol 8 stearate type I	11K04-N06	1
Noweda	Macrogol 8 stearate type I	12125404	1
Caelo	Macrogol 8 stearate type I	12125404	1
Caelo	Macrogol 8 stearate type I	12E24-N07	1
Fagron	Macrogol 8 stearate type I	12E24-N07	1
Caelo	Macrogol 8 stearate type I	12H23-N02	1
Caelo	Macrogol 8 stearate type I	13077505	2
Caelo	Macrogol 8 stearate type I	13077501	1
Fagron	Macrogol 8 stearate type I	T30911	1
Fagron	Macrogol 8 stearate type I	1115M-03130	1
Fagron	Macrogol 8 stearate type I	14H04-B01-298920	1

- 91 303 spectra from 501 *Apo-Ident* customers from a total of 26 985 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Macrogol 8 stearate type I* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Macrogol 8 stearate type I* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	120	0	93 754
Type C	0	48	1	91 303

The substance/substance group *Macrogol 8 stearate type I* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6587 %)	97.9592 % (> 91.8367 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31165	31165	0.00	20.42
31453	31453	0.00	28.34
33717	33717	0.00	39.13
33900	33900	0.00	20.18

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Macrogol ointment**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30465-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Macrogol ointment; Macrogoli unguentum DAC; Unguentum polyäthylenglycoli

Special notes

When selecting the *Macrogol ointment* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Macrogol ointment*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Macrogol ointment	13002502	31004	40	20130108
Caelo	Macrogol ointment	13177801	31247	40	1405327
Caelo	Macrogol ointment	14078601	31442	40	20140321
Caelo	Macrogol ointment	13177801	31486	40	20130709

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Macrogol ointment*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 384 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 300 spectra from 7 reference samples of the substance/substance group *Macrogol ointment*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Macrogol ointment	10033301	30255	40
Caelo	Macrogol ointment	12005903	30465	40
Caelo	Macrogol ointment	153838	31989	60
Caelo	Macrogol ointment	162577	32493	40
Caelo	Macrogol ointment	162578	32494	40
Caelo	Macrogol ointment	172237	33751	40
Caelo	Macrogol ointment	182006	34257	40

- 93 574 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 147 spectra from 65 *Apo-Ident* customers from 52 batches from the substance/substance group *Macrogol ointment*.
- These include spectra of independent samples from 48 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macrogol ointment	6031507	1
Caelo	Macrogol ointment	12179102	4
Caelo	Macrogol ointment	1012A-02834	1
Fagron	Macrogol ointment	12661402	1
Caelo	Macrogol ointment	13177801	7
Caelo	Macrogol ointment	13324902	1
Caelo	Macrogol ointment	14012201	6
Fagron	Macrogol ointment	14012201	1
Caelo	Macrogol ointment	14078601	9
Caelo	Macrogol ointment	14219802	3
Caesar & Loretz GmbH	Macrogol ointment	14078601	1
Caelo	Macrogol ointment	14219801	1
Caelo	Macrogol ointment	14336301	2
Caesar & Loretz GmbH	Macrogol ointment	14336301	1
Caelo	Macrogol ointment	14336302	3
Caelo	Macrogol ointment	15170602	2
Caelo	Macrogol ointment	15170701	4
Caelo	Macrogol ointment	151707107	1
Caelo	Macrogol ointment	15383802	1
Caelo	Macrogol ointment	15383902	2
GEHE	Macrogol ointment	16186301	1
Caelo	Macrogol ointment	16186301	9
Caelo	Macrogol ointment	16062002	2
Caelo	Macrogol ointment	16062003	4
Caelo	Macrogol ointment	16092003	1
Caelo	Macrogol ointment	16186401	1
Caelo	Macrogol ointment	16257801	3
Caelo	Macrogol ointment	16257802	2
Caelo	Macrogol ointment	16314502	5
GEHE	Macrogol ointment	16314502	1
Caelo	Macrogol ointment	17127601	5
Caelo	Macrogol ointment	17063901	1
Caelo	Macrogol ointment	4511M-02834	1
Caelo	Macrogol ointment	6091202	1
GEHE	Macrogol ointment	17127601	2
Caelo	Macrogol ointment	3315A-02834	1
Caelo	Macrogol ointment	312E-02834	1
Ebert+Jacobi	Macrogol ointment	11277701	1
Caelo	Macrogol ointment	11277702	1
cesar&lorenz/Phönix	Macrogol ointment	11327702	1
Caelo	Macrogol ointment	12117202	2
Großhandel	Macrogol ointment	12179101	2
Anzag/	Macrogol ointment	12261402	2
Noweda	Macrogol ointment	12179102	1
Caelo	Macrogol ointment	12261402	1
Anzag	Macrogol ointment	12261402	1
Caelo	Macrogol ointment	12353501	5
Caelo	Macrogol ointment	12261403	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogol ointment	12353502	4
Anzag	Macrogol ointment	12353501	3
Caelo	Macrogol ointment	13002502	7
Caelo	Macrogol ointment	13002602	5
Großhandel	Macrogol ointment	13002602	1
Caelo	Macrogol ointment	13177701	1
Caelo	Macrogol ointment	13177702	3
Bombastus	Macrogol ointment	13177801	1
Sanacorp	Macrogol ointment	13177801	1
Caelo	Macrogol ointment	13275302	1
Caelo	Macrogol ointment	13324901	1
Caelo	Macrogol ointment	22051304	1
Caelo	Macrogol ointment	19071305	1
Caelo	Macrogol ointment	29111202	1
Caelo	Macrogol ointment	10166302	1
Euro OTC	Macrogol ointment		1
Anzag	Macrogol ointment	12117202	1
Caelo	Macrogol ointment	13002601	1

- 91 205 spectra from 501 *Apo-Ident* customers from a total of 26 959 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Macrogol ointment* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Macrogol ointment* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	300	0	93 574
Type C	0	137	10	91 205

The substance/substance group *Macrogol ointment* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.0000 %)
Type C	100.0000 % (> 99.6586 %)	93.1973 % (> 91.1565 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31004	31247	19.83	38.02
31247	31247	0.00	42.76
31442	31247	17.71	34.02
31486	31247	9.90	37.06

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Macrogol-2-laurylether**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31731-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Macrogol-2-laurylether; Dehydol® LS 2

Special notes

When selecting the *Macrogol-2-laurylether* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Macrogol-2-laurylether*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Macrogol-2-laury...	14054304	31731	40	1603434
Caelo	Macrogol-2-laury...	14054309	31867	40	20140225

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Macrogol-2-laurylether*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 104 spectra from 4 reference samples of the substance/substance group *Macrogol-2-laurylether*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Macrogol-2-laurylether	14054304	31731 [†]	20
Caelo	Macrogol-2-laurylether	14054309	31867 [†]	20
Caelo	Macrogol-2-laurylether	15411908	33273	40
Caelo	Macrogol-2-laurylether	15411904	not required	24

- 93 770 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 4 spectra from 4 *Apo-Ident* customers from 3 batches from the substance/substance group *Macrogol-2-laurylether*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macrogol-2-laurylether	15411903	1
Caelo	Macrogol-2-laurylether	15411904	1
Caelo	Macrogol-2-laurylether	15411902	2

- 91 348 spectra from 501 *Apo-Ident* customers from a total of 27 006 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Macrogol-2-laurylether* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Macrogol-2-laurylether* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	104	0	93 770
Type C	0	4	0	91 348

The substance/substance group *Macrogol-2-laurylether* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.2308 %)
Type C	100.0000 % (> 99.6608 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31731	31731	0.00	25.58
31867	31731	7.08	23.25

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Macrogols 300**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31151-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Macrogols 300; Polyethylenglycol 300

Special notes

When selecting the *Macrogols 300* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Macrogols 300*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Macrogols 300	13070109	31151	40	1403339
Caelo	Macrogols 300	14043802	31370	40	20140221
Caelo	Macrogols 300	16022102	32073	40	20160222

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Macrogols 300*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 260 spectra from 7 reference samples of the substance/substance group *Macrogols 300*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Macrogols 300	09375703	30152	40
Caelo	Macrogols 300	16022102	32073 [†]	20
Caelo	Macrogols 300	160221	32209	40
Caelo	Macrogols 300	161609	32212	40
Caelo	Macrogols 300	162196	32337	40
Caelo	Macrogols 300	172105	33710	40
Caelo	Macrogols 300	181748	34232	40

- 93 614 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 40 spectra from 23 *Apo-Ident* customers from 28 batches from the substance/substance group *Macrogols 300*.
- These include spectra of independent samples from 25 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 300	1323304	1
Caelo	Macrogols 300	13233304	5
Phoenix am 13.03.14	Macrogols 300	13233304	1
Caelo	Macrogols 300	14145106	2
Caelo	Macrogols 300	14043802	2
Caelo	Macrogols 300	14145107	1
Caelo	Macrogols 300	14183102	2
Caelo	Macrogols 300	15286402	2
Caelo	Macrogols 300	17034401	1
Caelo	Macrogols 300	17034406	1
Euro OTC	Macrogols 300	4315A-02102	1
Euro OTC	Macrogols 300	5214I-02102	1
Caelo	Macrogols 300	11043206	1
Anzag 10.12.2012	Macrogols 300	11170206	1
Caelo	Macrogols 300	11336204	1
Caelo	Macrogols 300	11336201	2
Caelo	Macrogols 300	11336205	1
Caelo	Macrogols 300	1133624	1
Caelo	Macrogols 300	12137803	1
Caelo	Macrogols 300	12108101	1
Gehe	Macrogols 300	12137805	1
Caelo	Macrogols 300	12108103	1
Caelo	Macrogols 300	12137805	1
Caelo	Macrogols 300	12272906	1
Caelo	Macrogols 300	1307109	1
Gehe 21.06.2013	Macrogols 300	13070109	1
Caelo	Macrogols 300	13070109	1
Caelo	Macrogols 300	1406212102	1
Euro OTC	Macrogols 300	2512051	1
Caelo	Macrogols 300		1
Euro OTC	Macrogols 300	4312I-02102	1

- 91 312 spectra from 501 *Apo-Ident* customers from a total of 26 981 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Macrogols 300* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Macrogols 300* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following

table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	260	0	93 614
Type C	0	37	3	91 312

The substance/substance group *Macrogols 300* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.6923 %)
Type C	100.0000 % (> 99.6588 %)	92.5000 % (> 85.0000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31151	31151	0.00	51.12
31370	31151	4.11	47.33
32073	31151	12.77	43.25

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Macrogols 400**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30523-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Macrogols 400; Macrogolum 400; Polyäthylenglycol 400

Special notes

When selecting the *Macrogols 400* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Macrogols 400*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Macrogols 400	13171112	31149	40	1403644
Caelo	Macrogols 400	14009101	31367	40	20140123
Euro OTC	Macrogols 400	1510042	31814	40	20151216*
Euro OTC	Macrogols 400	1609017	32714	40	20160929*
Caelo	Macrogols 400	170137	33057	40	20170125*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Macrogols 400*. These samples are listed above in the section *calibration samples*. The reference samples originate from 5 different batches.
- 35 344 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 580 spectra from 15 reference samples of the substance/substance group *Macrogols 400*.
- These include spectra of independent samples from 14 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Macrogols 400	10105810	30166	40
Euro OTC	Macrogols 400	1008004-02	30523	40
Caelo	Macrogols 400	12043301	30524	40
Euro OTC	Macrogols 400	1510042	31814 [†]	20
Euro OTC	Macrogols 400	1603008	32185	40
Euro OTC	Macrogols 400	1412015	32187	40
Caelo	Macrogols 400	160022	32211	40
Euro OTC	Macrogols 400	1607022	32461	40
Euro OTC	Macrogols 400	1610001	32712	40
Caelo	Macrogols 400	162941	33096	40
Euro OTC	Macrogols 400	1702024	33226	40
Caelo	Macrogols 400	171100	33321	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Macrogols 400	172833	33890	40
Euro OTC	Macrogols 400	1711029	33932	40
Euro OTC	Macrogols 400	1807027	34214	40

- 93 294 spectra from a total of 2152 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 505 spectra from 176 *Apo-Ident* customers from 145 batches from the substance/substance group *Macrogols 400*.
- These include spectra of independent samples from 143 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 400	12152602	4
Caelo	Macrogols 400	12152614	3
Caelo	Macrogols 400	3790505	1
Fafron/Anzag	Macrogols 400	1014E-03059	1
Caelo	Macrogols 400	12277607	9
Fagron	Macrogols 400	13/26-N21	1
Euro OTC	Macrogols 400	130403001	1
Euro OTC	Macrogols 400	1304030-01	10
Caelo	Macrogols 400	13171103	5
Fagron	Macrogols 400	13126-N21	1
AHD	Macrogols 400	13171106	1
Caelo	Macrogols 400	13171107	5
Caelo	Macrogols 400	13171112	4
Caelo	Macrogols 400	13272901	5
Caelo	Macrogols 400	13272906	4
Caelo	Macrogols 400	13272904	6
Alliance Health 14.04.2...	Macrogols 400	13272908	2
Caelo	Macrogols 400	13272909	3
Caelo	Macrogols 400	13272908	1
Fagron	Macrogols 400	13F06-B40-289978	1
Fagron	Macrogols 400	13F06-N11	10
Caelo	Macrogols 400	14009101	3
Caelo	Macrogols 400	14009102	4
Caelo	Macrogols 400	14009103	4
Caelo	Macrogols 400	14009107	6
Caelo	Macrogols 400	14009109	5
Phönix	Macrogols 400	14009109	1
Caelo	Macrogols 400	14009110	4

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 400	14009113	3
Euro OTC	Macrogols 400	1402011-01	13
Euro OTC	Macrogols 400	1402011-11	1
Caelo	Macrogols 400	1409103	1
Caelo	Macrogols 400	140211-01	1
OTG/Phönix	Macrogols 400	1412015-01	1
Euro OTC	Macrogols 400	1412015-01	5
Caelo	Macrogols 400	1412015-01	1
Caelo	Macrogols 400	14195601	18
Euro OTC	Macrogols 400	1412015-02	1
Caelo	Macrogols 400	14195602	8
Caelo	Macrogols 400	14195603	12
Caesar & Loretz GmbH	Macrogols 400	14195603	1
Caelo	Macrogols 400	14195605	5
Caelo	Macrogols 400	14195613	1
Caelo	Macrogols 400	14352801	4
Caelo	Macrogols 400	14352803	13
Caesar&Loretz GmbH/Anzag	Macrogols 400	14352803	2
Caelo	Macrogols 400	14352804	4
Canea/Phönix	Macrogols 400	14352805	1
Caelo	Macrogols 400	14352805	2
Caelo	Macrogols 400	14352811	1
Caelo	Macrogols 400	14352813	3
Caelo	Macrogols 400	14352814	14
Bombastus	Macrogols 400	14352814	2
Phönix	Macrogols 400	14352814	1
Fagron	Macrogols 400	14B05-B04	1
Fagron	Macrogols 400	14B05-B04-294688	3
Fagron	Macrogols 400	14C10-B05-300393	1
Fagron	Macrogols 400	14G10-B05	2
Fagron	Macrogols 400	14G10-B05300393	1
Fagron	Macrogols 400	14G10-B05-300393	2
Fagron	Macrogols 400	14G10-B05-300394	2
Fagron	Macrogols 400	14K26-B01	2
Fagron	Macrogols 400	14K26-B01-301176	1
Euro OTC	Macrogols 400	1510042-01	4
Caelo	Macrogols 400	15131401	2
Caelo	Macrogols 400	15131406	2
Noweda	Macrogols 400	15131406	1
Hustadt-Apotheke	Macrogols 400	15131407	1
Caelo	Macrogols 400	15324001	5
Caelo	Macrogols 400	15324002	15
Caelo	Macrogols 400	15324003	1
Caelo	Macrogols 400	15324004	1
Hustadt-Apotheke	Macrogols 400	15324003	1
Caelo	Macrogols 400	15324007	1
Noweda	Macrogols 400	15324008	1
Caelo	Macrogols 400	15324008	2
Caelo	Macrogols 400	16002201	20
Fagron	Macrogols 400	15B16-B02-308528	1
Caelo	Macrogols 400	15H28-B11-313603	1
Gehe	Macrogols 400	16002201	1
Caelo	Macrogols 400	16002202	1
Caelo	Macrogols 400	16002203	1
Phoenix	Macrogols 400	16002203	1
Euro OTC	Macrogols 400	1603008-01	3
Euro OTC	Macrogols 400	1607022-01	2
Euro OTC	Macrogols 400	1610001-01	1
Bombastus	Macrogols 400	160300801	1
Caelo	Macrogols 400	16128203	19

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
WE:Sanacorp 19.10.2016 ...	Macrogols 400	16128203	1
Caelo	Macrogols 400	16128204/170227M	1
Caelo	Macrogols 400	16128205	6
phönix	Macrogols 400	16128205	1
Caelo	Macrogols 400	16294102	5
Caelo	Macrogols 400	16294106	3
Caelo	Macrogols 400	16294107	1
Finteler Apotheke	Macrogols 400	17013701	1
Caelo	Macrogols 400	17013701	3
Caelo	Macrogols 400	17013702	2
Caelo	Macrogols 400	17013706	4
Euro OTC	Macrogols 400	1702024-01	1
Euro OTC	Macrogols 400	1702024-02	1
Caelo	Macrogols 400	17110001	4
Fagron	Macrogols 400	17110001	1
Finteler Apotheke	Macrogols 400	17110001	1
Caelo	Macrogols 400	17110007	1
Caelo	Macrogols 400	17110006	1
Caelo	Macrogols 400	17283301	1
Euro OTC	Macrogols 400	17112019-01	1
Hustadt-Apotheke	Macrogols 400	1716I-03059	1
Caelo	Macrogols 400	17283302	2
Euro OTC	Macrogols 400	180308M/1711029-01	1
Caelo	Macrogols 400	18081505	1
Caelo	Macrogols 400	21092015B	1
Fagron	Macrogols 400	2314A-03059	1
Caelo	Macrogols 400	2412U-03059	1
Euro OTC	Macrogols 400	1008004-02	1
Fagron	Macrogols 400	11/45	1
Caelo	Macrogols 400	11293302	1
Bombastus	Macrogols 400	3215E-03059	1
Caelo	Macrogols 400	11172707	1
Caelo	Macrogols 400	11172705	1
Caelo	Macrogols 400	11293304	1
Caelo	Macrogols 400	11293311	1
Noweda	Macrogols 400	11293312	1
Fagron	Macrogols 400	11F07-N01	1
Fagron	Macrogols 400	11I13-N04	1
Caelo	Macrogols 400	12043301	3
Cealo/Gehe	Macrogols 400	12043311	1
Caelo	Macrogols 400	12043311	1
Phönix 06.10.2012	Macrogols 400	12043306	1
Caelo	Macrogols 400	12043312	2
Noweda	Macrogols 400	12043312	1
Euro OTC	Macrogols 400	1206030-01	3
Euro OTC	Macrogols 400	1206030-02	5
Phönix	Macrogols 400	1206030-01	1
Caelo	Macrogols 400	1213M-03059	1
Caelo	Macrogols 400	12152601	1
Anzag	Macrogols 400	12152601	1
Caelo	Macrogols 400	12152608	7
Caelo	Macrogols 400	12152607	2
Hustadt Apotheke	Macrogols 400	12152607	1
Caelo	Macrogols 400	12152612	8
Caelo	Macrogols 400	12152611	3
Phoenix	Macrogols 400	12152611	1
Caelo	Macrogols 400	12152613	4
Noweda	Macrogols 400	12152614	3
Phönox 26.04.2013	Macrogols 400	12152614	1
Noweda	Macrogols 400	1215614	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 400	12277603	2
Phönix Gotha	Macrogols 400	12277607	1
Caelo	Macrogols 400	12277608	1
Caelo	Macrogols 400	12277611	1
Gehe 25.04.2013	Macrogols 400	12293312	1
Fagron	Macrogols 400	12C21-N10	7
Fagron	Macrogols 400	12c21-n10	1
Fagron	Macrogols 400	12H02-N05	5
Euro OTC	Macrogols 400	1304030	1
Fagron	Macrogols 400	12H02-N05	1
Caelo	Macrogols 400	13171102	9
Caelo	Macrogols 400	13171104	5
Caelo	Macrogols 400	13171106	4
Caesar/Noweda	Macrogols 400	13171106	1
Caelo	Macrogols 400	131711102	1
Caelo	Macrogols 400	13171109	3
Caelo	Macrogols 400	13171113	2
Fagron	Macrogols 400	13A25-N02	7
Fagron	Macrogols 400	13I26-N21	1
Hedinger	Macrogols 400	13F06-N11	1
Caelo	Macrogols 400	27111309	1
Caelo	Macrogols 400	20021306	1
Caelo	Macrogols 400	20190612-1	1
Caelo	Macrogols 400	38121-03059	1
AHD;19.06.15;3,34EUR	Macrogols 400	14352801	1
sanacorp, 3,43EUR, 20.09...	Macrogols 400	16128201	1
Caelo	Macrogols 400	1314I-03059	1
VDL 16.04.14	Macrogols 400	13272906//05R06	1
Caelo	Macrogols 400	1206030-02	2
Phönix	Macrogols 400	12152601	1

- 90 847 spectra from 501 *Apo-Ident* customers from a total of 26 867 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Macrogols 400* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Macrogols 400* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	570	10	93 294
Type C	0	487	18	90 847

The substance/substance group *Macrogols 400* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	98.2759 % (> 97.7586 %)
Type C	100.0000 % (> 99.6586 %)	96.4356 % (> 95.8416 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31149	31149	0.00	23.95
31367	33057	8.18	20.46
31814	31814	0.00	24.93
32714	32714	0.00	22.39
33057	33057	0.00	23.14

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Mandarin oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31747-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Mandarin oil; Oleum Mandarinae

Special notes

When selecting the *Mandarin oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Mandarin oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Mandarin oil	14064103	31747	40	1602210
Caelo	Mandarin oil	15398001	31863	40	20151127

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Mandarin oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 132 spectra from 5 reference samples of the substance/substance group *Mandarin oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Mandarin oil	14064103	31747 [†]	20
Caelo	Mandarin oil	15398001	31863 [†]	20
Caelo	Mandarin oil	170368	33074	40
Caelo	Mandarin oil	180923	34055	40
Caelo	Mandarin oil	15011601	not required	12

- 93 742 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 5 spectra from 2 *Apo-Ident* customers from 4 batches from the substance/substance group *Mandarin oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Mandarin oil	15398001	1
Caelo	Mandarin oil	15398003	1
Caelo	Mandarin oil	17036801	2
Taoasis	Mandarin oil	412011-124957	1

- 91 347 spectra from 501 *Apo-Ident* customers from a total of 27 005 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Mandarin oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Mandarin oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	132	0	93 742
Type C	0	5	0	91 347

The substance/substance group *Mandarin oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.4545 %)
Type C	100.0000 % (> 99.6603 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31747	31747	0.00	19.97
31863	31747	5.18	19.18

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Matricaria oil, rich in (-)-alpha-bisabolol
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30738-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Matricaria oil, rich in (-)-alpha-bisabolol; Oleum chamomillae

Special notes

When selecting the *Matricaria oil, rich in (-)-alpha-bisabolol* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Matricaria oil, rich in (-)-alpha-bisabolol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Matricaria oil, ...	11221206	30738	40	20110810
Caelo	Matricaria oil, ...	13195103	31239	30	1504133
Caelo	Matricaria oil, ...	13334702	31427	40	20131009
Caelo	Matricaria oil, ...	13334702	31428	40	20131009

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 150 spectra of 4 reference samples from the substance/substance group *Matricaria oil, rich in (-)-alpha-bisabolol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 394 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 5 reference samples of the substance/substance group *Matricaria oil, rich in (-)-alpha-bisabolol*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Matricaria oil, rich in (-)-...	13195101	31160	20
Caelo	Matricaria oil, rich in (-)-...	13334702	31427 [†]	20
Caelo	Matricaria oil, rich in (-)-...	13334702	31428 [†]	20
Caelo	Matricaria oil, rich in (-)-...	163310	33081	40
Caelo	Matricaria oil, rich in (-)-...	181173	34072	40

- 93 734 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 12 spectra from 6 *Apo-Ident* customers from 10 batches from the substance/substance group *Matricaria oil, rich in (-)-alpha-bisabolol*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Matricaria oil, rich in (-)-...	15080505	2
Caelo	Matricaria oil, rich in (-)-...	15137801	2
Caelo	Matricaria oil, rich in (-)-...	18091401	1
Taoasis	Matricaria oil, rich in (-)-...	16266-119132	1
Taoasis	Matricaria oil, rich in (-)-...	6181024-122141	1
Taoasis	Matricaria oil, rich in (-)-...	6181024-121466	1
Taoasis	Matricaria oil, rich in (-)-...	6181024-123770	1
Taoasis	Matricaria oil, rich in (-)-...	929-115334	1
Caelo	Matricaria oil, rich in (-)-...	16048910	1
Taoasis	Matricaria oil, rich in (-)-...	6181024-123496	1

- 91 340 spectra from 501 *Apo-Ident* customers from a total of 26 999 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Matricaria oil, rich in (-)-alpha-bisabolol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Matricaria oil, rich in (-)-alpha-bisabolol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	150	0	35 394
Type B	0	140	0	93 734
Type C	0	12	0	91 340

The substance/substance group *Matricaria oil, rich in (-)-alpha-bisabolol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6593 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30738	31239	15.00	42.95
31239	31239	0.00	42.51
31427	31239	6.71	41.42
31428	31239	6.90	41.48

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Melissa oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30207-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Melissa oil; Oleum melissae rectific.

Special notes

When selecting the *Melissa oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	European Pharmacopoeia 8th Edition, Basic Version 2014 [3]
Komm2.2.40	Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident [4]
AA004	Erstellung und Validierung eines IdentModul-Updates
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Melissa oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Melissa oil	13086306	31163	40	1403319
Caelo	Melissa oil	13354110	31403	40	20131028
Caelo	Melissa oil	163349	33080	40	20161215*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Melissa oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 130 spectra from 3 reference samples of the substance/substance group *Melissa oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Melissa oil	10019701	30207	50
Caelo	Melissa oil	16053609	33784	40
Caelo	Melissa oil	180973	34084	40

- 93 744 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 131 spectra from 62 *Apo-Ident* customers from 48 batches from the substance/substance group *Melissa oil*.
- These include spectra of independent samples from 46 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Melissa oil	12283804	6
Caelo	Melissa oil	12283809	8
Caelo	Melissa oil	12253804	1
Nachprüfung	Melissa oil	130863060	1
Caelo	Melissa oil	13086308	6
Caelo	Melissa oil	13086311	6
Caelo	Melissa oil	13354110	3
Sanacorp	Melissa oil	1334110	1
Caelo	Melissa oil	13354113	5
Caelo	Melissa oil	14332511	13
Caelo	Melissa oil	14332501	5
Caelo	Melissa oil	1512303	1
Caelo	Melissa oil	15172303	9
Caelo	Melissa oil	15172308	6
Caelo	Melissa oil	18081501	1
Caelo	Melissa oil	11035803	1
Caelo	Melissa oil	11323902	3
Caesar&lorenz,Anzag	Melissa oil	11323902	1
Caesar & Lorenz GmbH, 4...	Melissa oil	11323910	1
Caelo	Melissa oil	11323903	1
Caelo	Melissa oil	122283804	1
Caelo	Melissa oil	11323911	1
Anzag 27.02.13 ek: 8.17EUR	Melissa oil	122833801	1
Caelo	Melissa oil	12283801	2
celo/phönix	Melissa oil	12283805	1
Caelo	Melissa oil	12283805	2
Phönix	Melissa oil	12283809	1
Gehe 14.06.2013	Melissa oil	12283809	2
Caelo	Melissa oil	13086306	4
Sanacorp	Melissa oil	13086306	1
Caelo	Melissa oil	14332508	2
Caelo	Melissa oil	14332502	1
Caelo	Melissa oil	16053601	5
Alte Apotheke	Melissa oil	16053601	1
Caleo/Phönix	Melissa oil	16053609	1
Caelo	Melissa oil	16334905	1
Caelo	Melissa oil	16053609	2
Caelo	Melissa oil	17012014	1
Caelo	Melissa oil	17091201	1
Caelo	Melissa oil	1712070M/16053609	1
Caelo	Melissa oil	190412	1
Caelo	Melissa oil	20230712-7	1
Caelo	Melissa oil	21111303	1
Caelo	Melissa oil	5211R-01993	1
Caelo	Melissa oil	Ch12283804	1
Caelo	Melissa oil	28051311	1
Caelo	Melissa oil	13354105	2
Caelo	Melissa oil	13086310	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Melissa oil	14332504	2
Caelo	Melissa oil	13354108	1
Caelo	Melissa oil	14332512	1
Caelo	Melissa oil	13086304	1
Caelo	Melissa oil	15172310	1
Caelo	Melissa oil	4110I-01993	1
Caelo	Melissa oil	11323909	2

- 91 221 spectra from 501 *Apo-Ident* customers from a total of 26 962 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Melissa oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Melissa oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	130	0	93 744
Type C	0	129	2	91 221

The substance/substance group *Melissa oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.3846 %)
Type C	100.0000 % (> 99.6586 %)	98.4733 % (> 96.1832 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31163	31163	0.00	47.04
31403	31163	4.84	44.64
33080	33080	0.00	48.42

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Methyl salicylate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31038-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Methyl salicylate; Methylis salicylas; Methylum salicylicum

Special notes

When selecting the *Methyl salicylate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Methyl salicylate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Methyl salicylate	12215217	31038	40	1403162
Fagron	Methyl salicylate	14D03-B04-294169	31607	40	20140414
Fagron	Methyl salicylate	16E23-B05-332971	33275	40	20160608

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Methyl salicylate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 141 spectra from 4 reference samples of the substance/substance group *Methyl salicylate*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Methyl salicylate	10179004	30330	41
Fagron	Methyl salicylate	14D03-B04-294169	31607 [†]	20
Caelo	Methyl salicylate	15409810	33280	40
Fagron	Methyl salicylate	18B19-B06-346255	34047	40

- 93 733 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 78 spectra from 39 *Apo-Ident* customers from 47 batches from the substance/substance group *Methyl salicylate*.
- These include spectra of independent samples from 46 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Methyl salicylate	13185101	3
Bomb/ Noweda	Methyl salicylate	RMS1332311	1
Caelo	Methyl salicylate	12215215	1
Fagron	Methyl salicylate	13I11-N06	1
Caelo	Methyl salicylate	12308307	1
Caelo	Methyl salicylate	13189309	2
Caelo	Methyl salicylate	15117910	1
Fagron	Methyl salicylate	15809-B04	1
Fagron	Methyl salicylate	15B09-B04	1
Fagron	Methyl salicylate	15H04-B01	1
Fagron	Methyl salicylate	17H21-B04-338848	1
Caelo	Methyl salicylate	10313708	3
Fagron	Methyl salicylate	11H30-N14	1
Nachprüfung	Methyl salicylate	12D27-N05	2
Fagron	Methyl salicylate	1139	1
Fagron	Methyl salicylate	13/11-N06	4
Caelo	Methyl salicylate	12308313	1
Alliance Health 18.11.2...	Methyl salicylate	13185102	1
Caelo	Methyl salicylate	13189308	1
PHOENIX 08.05.2014	Methyl salicylate	13I11-N06	1
Fagron	Methyl salicylate	14D03-B04	1
Fagron	Methyl salicylate	14D03-B04-292643	1
Fagron	Methyl salicylate	14D03-B04-294169	2
Caelo	Methyl salicylate	14D03-B04-294169	1
Fagron	Methyl salicylate	14D03-B04-294170	1
Caelo	Methyl salicylate	15117904	2
Caelo	Methyl salicylate	15409807	1
Fagron	Methyl salicylate	15B09-04-309068	1
Caelo	Methyl salicylate	15B09-B04	1
Fagron	Methyl salicylate	15B09-B04303941	1
Fagron	Methyl salicylate	16e23-b05-332971	1
Caelo	Methyl salicylate	17071301	1
Fagron	Methyl salicylate	30121414	1
Fagron	Methyl salicylate	303941	1
ACM	Methyl salicylate	KNC1794J14	1
Fagron	Methyl salicylate	RMS1332311	1
Fagron	Methyl salicylate	RMS13332311	1
Alliance Health 10.06.2...	Methyl salicylate	RMS1332311	1
Fagron	Methyl salicylate	RMS1409712	1
Caelo	Methyl salicylate	7111203	1
Fagron	Methyl salicylate	711Q-01867	1
Fagron	Methyl salicylate	11D08-N07	3
Fagron	Methyl salicylate	12D24-N09	1
Fagron	Methyl salicylate	12C02-No2	1
Fagron	Methyl salicylate	12D27-n05	1
Fagron	Methyl salicylate	12D27-N05	3
Caesar & Loretz GmbH, 4...	Methyl salicylate	12215215	1
Fagron	Methyl salicylate	13/11-N060013/364	1
Fagron	Methyl salicylate	13/11-N060025/364	1
Fagron	Methyl salicylate	13D24-N09	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Methyl salicylate	13D24-N09	1
Fagron	Methyl salicylate	13d24-no9	1
Caelo	Methyl salicylate	13I11-N06	1
Fagron	Methyl salicylate	19061309	1
Fagron	Methyl salicylate	172650-0002	2
Caelo	Methyl salicylate	90328419	1

- 91 274 spectra from 501 *Apo-Ident* customers from a total of 26 963 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Methyl salicylate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Methyl salicylate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	141	0	93 733
Type C	0	78	0	91 274

The substance/substance group *Methyl salicylate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7447 %)
Type C	100.0000 % (> 99.6587 %)	100.0000 % (> 92.3077 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31038	31038	0.00	145.44
31607	31038	6.74	146.08
33275	31038	2.37	147.08

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Metronidazole 25% in white petrolatum**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30663-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Metronidazole 25% in white petrolatum

Special notes

When selecting the *Metronidazole 25% in white petrolatum* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Metronidazole 25% in white petrolatum*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Metronidazole 25...	12266401	30663	19	1405100
Caelo	Metronidazole 25...	14116801	31316	20	20140411
Caelo	Metronidazole 25...	16063101	32456	40	20160229
Caelo	Metronidazole 25...	16063101	32768	40	20160229

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 119 spectra of 4 reference samples from the substance/substance group *Metronidazole 25% in white petrolatum*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 425 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 4 reference samples of the substance/substance group *Metronidazole 25% in white petrolatum*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Metronidazole 25% in white p...	11214101	30420	40
Caelo	Metronidazole 25% in white p...	12266401	30663 [†]	20
Caelo	Metronidazole 25% in white p...	14116801	31316 [†]	20
Caelo	Metronidazole 25% in white p...	16249701	33272	40

- 93 754 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 39 spectra from 12 *Apo-Ident* customers from 14 batches from the substance/substance group *Metronidazole 25% in white petrolatum*.
- These include spectra of independent samples from 12 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Metronidazole 25% in white p...	14024301	4
Caelo	Metronidazole 25% in white p...	14116801	3
Caelo	Metronidazole 25% in white p...	14209301	1
Gehe	Metronidazole 25% in white p...	14209301	1
Caelo	Metronidazole 25% in white p...	14247201	1
Noweda 20.02.2015	Metronidazole 25% in white p...	14247301	1
Caelo	Metronidazole 25% in white p...	15207101	5
Cael Sanacorp 05.02.2016	Metronidazole 25% in white p...	15207101	1
Caelo	Metronidazole 25% in white p...	16063101	4
Caelo	Metronidazole 25% in white p...	13109801	4
celo/Jenne	Metronidazole 25% in white p...	13077601	1
Caelo	Metronidazole 25% in white p...	13188901	2
Caelo	Metronidazole 25% in white p...	13296501	3
Caelo	Metronidazole 25% in white p...	1285301	1
Sanacorp	Metronidazole 25% in white p...	12306601	1
GEHE/26.02.2013	Metronidazole 25% in white p...	12306601	1
Caelo	Metronidazole 25% in white p...	13002801	2
Caelo	Metronidazole 25% in white p...	13077601	3

- 91 313 spectra from 501 *Apo-Ident* customers from a total of 26 995 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Metronidazole 25% in white petrolatum* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Metronidazole 25% in white petrolatum* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	119	0	35 425
Type B	0	100	20	93 754
Type C	0	38	1	91 313

The substance/substance group *Metronidazole 25% in white petrolatum* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 94.9580 %)
Type B	100.0000 % (> 99.9761 %)	83.3333 % (> 80.8333 %)
Type C	100.0000 % (> 99.6588 %)	97.4359 % (> 89.7436 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30663	30663	0.00	55.97
31316	30663	6.92	51.09
32456	30663	16.08	56.75
32768	30663	23.30	57.76

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Mint oil / Peppermint oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30856-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Mint oil / Peppermint oil; Menthae arvensis aetheroleum partim mentholum depletum; Menthae piperitae aetheroleum; Mint oil; Oleum Menthae arvens. var. piper. japanisch; Oleum menthae arvensis; Oleum menthae piperitae dopp. rectific.; Peppermint oil

Special notes

When selecting the *Mint oil / Peppermint oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	European Pharmacopoeia 8th Edition, Basic Version 2014 [3]
Komm2.2.40	Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident [4]
AA004	Erstellung und Validierung eines IdentModul-Updates
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Mint oil / Peppermint oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Peppermint oil	12114206	30856	40	20120424

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Euro OTC	Peppermint oil	1312009-01	31294	40	20140130
Caelo	Peppermint oil	13303915	31410	40	1502136
Euro OTC	Peppermint oil	1512011	31811	40	20151214*
Caelo	Mint oil	16053503	32474	40	1701393
Caelo	Mint oil	15380507	32476	40	1701394
Caelo	Mint oil	16053503	32780	40	20160309
Caelo	Mint oil	15380509	32781	40	20151116

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 320 spectra of 8 reference samples from the substance/substance group *Mint oil / Peppermint oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 7 different batches.
- 35 224 spectra from a total of 833 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 300 spectra from 8 reference samples of the substance/substance group *Mint oil / Peppermint oil*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Peppermint oil	1512011	31811 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Peppermint oil	161904	32333	40
Euro OTC	Peppermint oil	1607018	32468	40
Caelo	Mint oil	15380507	32484	40
Caelo	Mint oil	16053508	32845	40
Euro OTC	Peppermint oil	1703032	33229	40
Caelo	Peppermint oil	172919	33699	40
Euro OTC	Peppermint oil	1804008	34089	40

- 93 574 spectra from a total of 2159 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 577 spectra from 250 *Apo-Ident* customers from 208 batches from the substance/substance group *Mint oil / Peppermint oil*.
- These include spectra of independent samples from 201 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Peppermint oil	13I04-N01	7
Bombastus	Peppermint oil	14150216	1
Fagron	Peppermint oil	13I04-N10	4
Fagron	Peppermint oil	13/04-N10	2
Fagron	Peppermint oil	14B03-B40-295665	3
Bombastus	Peppermint oil	285571	2
Bombastus	Peppermint oil	285715	2
Bombastus	Peppermint oil	286376	1
Bombastus	Peppermint oil	287325	1
Bombastus	Peppermint oil	287528	3
Bombastus	Peppermint oil	289675	1
Bombastus	Peppermint oil	290301	5
Caelo	Peppermint oil	1124709	1
Caelo	Peppermint oil	12114211	2
Fagron	Peppermint oil	700599	1
Fagron	Peppermint oil	12D10-N04	7
Fagron	Peppermint oil	13A10-N08	5
Bombastus	Peppermint oil	279930	5
Caelo	Peppermint oil	13303915	12
Caelo	Peppermint oil	13318004	4
Caelo	Peppermint oil	14150214	14
Caelo	Peppermint oil	14150221	4
Caelo	Peppermint oil	14369909	11
Sanacorp, 16.08.2016, 7...	Peppermint oil	15255816	2

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Peppermint oil	15255801	23
Caelo	Peppermint oil	15255811	8
Caelo	Peppermint oil	15255816	15
Fagron	Peppermint oil	15B11-B08-308854	1
Caelo	Peppermint oil	16114602	10
Bombastus	Peppermint oil	288864	2
Bombastus	Peppermint oil	292138	6
Bombastus	Peppermint oil	294293	3
Bombastus	Peppermint oil	295599	2
Bombastus	Peppermint oil	295737	4
Caelo	Peppermint oil	298969	1
Bombastus	Peppermint oil	299058	4
Caelo	Peppermint oil	11124706	1
Fagron	Peppermint oil	11g19-n07	1
Anzag	Peppermint oil	11G19N07	1
Fiebig	Peppermint oil	21022172/1	1
Bombastus	Peppermint oil	278812	2
Bombastus	Peppermint oil	274804	1
Bombastus	Peppermint oil	275481	3
Bombastus	Peppermint oil	280894	1
Bombastus	Peppermint oil	284588	1
Bombastus	Peppermint oil	286204	4
Caelo	Peppermint oil		2
Caelo	Peppermint oil	5170717	1
Primavera	Peppermint oil	5409142	1
Primavera	Peppermint oil	7703L22	1
Caelo	Peppermint oil	12180024	1
Caelo	Peppermint oil	13/04-N01	1
Caelo	Peppermint oil	11287407	6
Fagron	Peppermint oil	1304-N01	2
Caelo	Peppermint oil	13105801	10
Hedinger/Gehe	Peppermint oil	13105801	1
Gehe	Peppermint oil	13105801	1
Caelo	Peppermint oil	13241303	8
Caelo	Peppermint oil	13303903	5
Caelo	Peppermint oil	13303904	6
Caelo	Peppermint oil	13303912	2
Caelo	Peppermint oil	13303914	1
Sanacorp	Peppermint oil	13303915	3
Phönix	Peppermint oil	13303915	1
sanacorp	Peppermint oil	13I04-N01	1
Taoasis	Peppermint oil	13486-116309	1
Caelo	Peppermint oil	13I04-N01	1
Fagron	Peppermint oil	13104-N10	1
Euro OTC	Peppermint oil	1404038-01	1
Caeo	Peppermint oil	14150209	1
Euro OTC	Peppermint oil	1410005-01	1
Caelo	Peppermint oil	14150209	6
sanacorp	Peppermint oil	14150214	1
Caelo	Peppermint oil	14150215	2
Caelo	Peppermint oil	14150218	1
Caelo	Peppermint oil	14150220	6
Caelo	Peppermint oil	14150225	10
Phönix Gotha	Peppermint oil	14150225	1
Caelo	Peppermint oil	14369910	2
Caelo	Peppermint oil	14369911	2
Fagron	Peppermint oil	14B03-B40	3
Fagron	Peppermint oil	14B03B40	1
Euro OTC	Peppermint oil	1504017-01	3
Caelo	Peppermint oil	150415A	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Peppermint oil	15112908	1
Caelo	Peppermint oil	15112909	1
Euro OTC	Peppermint oil	1512011-01	5
Fagron	Peppermint oil	1512011-01	1
Caelo	Peppermint oil	15127004	3
Caelo	Peppermint oil	15126904	1
Caelo	Peppermint oil	15235504	2
Fagron	Peppermint oil	15255801	1
Caelo	Peppermint oil	15255805	2
Caesar/Noweda	Peppermint oil	15255811	1
Euro OTC	Peppermint oil	15255816	1
Caelo	Peppermint oil	15255817	4
Caelo	Peppermint oil	15363201	2
Caelo	Peppermint oil	15363207	1
Caelo	Peppermint oil	15363214	2
Caelo	Peppermint oil	1538707	1
Fagron	Peppermint oil	15B11-B08	1
Caelo	Peppermint oil	153663207	1
Caelo	Peppermint oil	15B11-B08-308854	1
Caelo	Peppermint oil	15383707	1
Fagron	Peppermint oil	15B11-B08-304117	1
Fagron	Peppermint oil	15K27-B01-317773	2
Euro OTC	Peppermint oil	1607018-01	1
Euro OTC	Peppermint oil	1607018-02	1
Caelo	Peppermint oil	16114601	1
Caelo	Peppermint oil	16114603	3
hochhaus-Apotheke	Peppermint oil	16114603	1
Caelo	Peppermint oil	16114604	1
Caelo	Peppermint oil	16114607	2
Caelo	Peppermint oil	16114611	8
Caelo	Peppermint oil	16114613	1
Sanacorp	Peppermint oil	16114625	1
Caelo	Peppermint oil	16114627	3
Caelo	Peppermint oil	16190405	4
Fagron	Peppermint oil	16F01-B03-322190	2
Alliance Healthcare	Peppermint oil	16F01-B03-322190	1
Fagron	Peppermint oil	16F01-B03-322818	1
Fagron	Peppermint oil	16F01-B03-326342	1
Euro OTC	Peppermint oil	1703032-01	1
Caelo	Peppermint oil	1718610	1
Caelo	Peppermint oil	17188609	2
Caelo	Peppermint oil	17188611	1
EuRho/Noweda	Peppermint oil	17291901	1
Fagron	Peppermint oil	17F08-B07-139332	1
GEHE	Peppermint oil	17F08-B07-324991	2
Fagron	Peppermint oil	17F08-B07-336043	1
Caelo	Peppermint oil	20111417	1
GEHE 06.06.2014	Peppermint oil	283491	1
Bombastus	Peppermint oil	283491	4
Bombastus	Peppermint oil	283734	1
Fagron	Peppermint oil	285571	1
Bombastus	Peppermint oil	289595	1
Bombastus	Peppermint oil	290394	2
Bombastus	Peppermint oil	292472	8
Bombastus	Peppermint oil	294802	1
Bombastus	Peppermint oil	295985	2
Bombastus	Peppermint oil	296017	1
Bombastus	Peppermint oil	297256	4
Caelo	Peppermint oil	297256	1
Bombastus	Peppermint oil	300468	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Peppermint oil	300843	2
Bombastus	Peppermint oil	302909	1
Bombastus	Peppermint oil	304515	4
Bombastus	Peppermint oil	303064	1
Bombastus	Peppermint oil	305627	1
Bombastus	Peppermint oil	309390	1
Taoasis	Peppermint oil	31145-118741	1
Bombastus	Peppermint oil	311547	1
Taoasis/Phönix	Peppermint oil	32724-118746	1
Taoasis	Peppermint oil	37979-120290	1
Caelo	Peppermint oil	3913I-01995	1
ACM	Peppermint oil	407000718901021	1
Caelo	Peppermint oil	4214U-01995	1
Caelo	Peppermint oil	52/AKD15	1
Fargon	Peppermint oil		1
Alliance	Peppermint oil	56335	1
Taoasis	Peppermint oil	6316-115825	1
Bio-Diät / Phönix	Peppermint oil	20553	1
Fargon	Peppermint oil	100030-N05	1
Bombastus	Peppermint oil	1001311995	1
Caelo	Peppermint oil	10235317	1
Caelo	Peppermint oil	11071312	1
Euro OTC	Peppermint oil	110800801	1
Caelo	Peppermint oil	11124709	3
Caelo	Peppermint oil	11124710	3
Noweda	Peppermint oil	11124710	1
Phönix	Peppermint oil	11124710	1
Caelo	Peppermint oil	11287416	9
Sanacorp	Peppermint oil	11287416	1
Caelo	Peppermint oil	11287419	1
Caelo	Peppermint oil	11287421	5
Fargon	Peppermint oil	11A05-N02	1
Fargon	Peppermint oil	11C14-N06	2
Fargon	Peppermint oil	11F09-N12	1
Fargon	Peppermint oil	11G19-N07	4
Fargon	Peppermint oil	12/19-N09	1
Caelo	Peppermint oil	12020803	1
Caelo	Peppermint oil	12020815	1
Caelo	Peppermint oil	12114206	15
Caelo	Peppermint oil	12020812	1
Fargon	Peppermint oil	12114206	1
Hedinger	Peppermint oil	12114206	1
Noweda	Peppermint oil	12114206	1
Sanacorp	Peppermint oil	12114206	1
Phoenix	Peppermint oil	12114211	1
Caelo	Peppermint oil	12114219	4
Sanacorp	Peppermint oil	12114219	1
Caelo	Peppermint oil	12114224	2
Caelo	Peppermint oil	12114225	15
Caesar & Loretz GmbH, 4...	Peppermint oil	12184601	2
Euro OTC	Peppermint oil	1303015-01	2
Caelo	Peppermint oil	13105306	1
Caelo	Peppermint oil	13105905	3
Caelo	Peppermint oil	1324302	1
Caelo	Peppermint oil	13241306	1
Gehe	Peppermint oil	13303903	1
Caesar/Noweda	Peppermint oil	13303915	1
Caelo	Peppermint oil	13318201	1
Fargon	Peppermint oil	13A10N08	2
Caelo	Peppermint oil	13A10-N08	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Peppermint oil	13A10-N8	1
Fagron	Peppermint oil	13i04-N10	1
Caelo	Peppermint oil	152555801	1
Caelo	Peppermint oil	20151112-6	1
Caelo	Peppermint oil	2112E-01995	1
Fagron	Peppermint oil	250214PJK02	1
Fagron	Peppermint oil	267849	1
Bombastus	Peppermint oil	269560	1
Anzag ek:13.70EUR 15.08.12	Peppermint oil	270566	1
Caelo	Peppermint oil	269560	1
Bombastus	Peppermint oil	271257	1
Phönix 10.01.13	Peppermint oil	273025	1
Bombastus	Peppermint oil	274976	1
Bombastus	Peppermint oil	277504	5
Bombastus	Peppermint oil	277817	2
DEFEKTUR	Peppermint oil	277504	1
Bombastus	Peppermint oil	282542	1
AHD, 29.11.13	Peppermint oil	282542	1
Sanacorp	Peppermint oil	280894	1
Anzag	Peppermint oil	283491	1
Phönix	Peppermint oil	287529	1
Bombastus	Peppermint oil	297946	1
phönix	Peppermint oil		1
Caelo	Peppermint oil	46/AD1	1
Sanacorp	Peppermint oil	47/AKD19	1
Taoasis	Peppermint oil	501319-105273	1
Taoasis	Peppermint oil	512777-107965	1
Taoasis	Peppermint oil	515574-107811	1
Taoasis	Peppermint oil	525728-109630	1
Taoasis GmbH	Peppermint oil	534640-110874	1
Taoasis	Peppermint oil	539490-111157	2
Taoasis	Peppermint oil	548992-112346	3
Caelo	Peppermint oil	688901	2
Fagron	Peppermint oil	688901	1
Bombastus	Peppermint oil	80312004	1
Fagron	Peppermint oil	8574769	2
Caelo	Peppermint oil	C101761	1
Taomed	Peppermint oil		1
Sensient	Peppermint oil	C9982	1
Caelo	Peppermint oil	Ch12114219	1
	Peppermint oil	9183299	1
Bombastus	Peppermint oil		1

- 90 775 spectra from 501 *Apo-Ident* customers from a total of 26 804 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Mint oil* / *Peppermint oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Mint oil* / *Peppermint oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or

incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	320	0	35 224
Type B	0	300	0	93 574
Type C	0	541	36	90 775

The substance/substance group *Mint oil / Peppermint oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.1250 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.0000 %)
Type C	100.0000 % (> 99.6586 %)	93.7608 % (> 93.2409 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30856	31410	11.52	54.53
31294	31811	2.62	55.08
31410	31410	0.00	56.96
31811	31811	0.00	57.03
32474	32474	0.00	56.91
32476	32476	0.00	54.25
32780	32474	3.80	56.80
32781	32476	5.40	55.50

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Myrtle oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31836-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Myrtle oil; Oleum myrte

Special notes

When selecting the *Myrtle oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Myrtle oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Myrtle oil	15225102	31836	40	1606483
Caelo	Myrtle oil	15225102	32349	40	20150721
Caelo	Myrtle oil	17227302	33879	40	1804344

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Myrtle oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 424 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Myrtle oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Myrtle oil	15225102	31836 [†]	20
Caelo	Myrtle oil	15225103	32182	40
Caelo	Myrtle oil	173184	34082	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2 spectra from 1 *Apo-Ident* customers from 2 batches from the substance/substance group *Myrtle oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Myrtle oil	15225104	1
Caelo	Myrtle oil	16296202	1

- 91 350 spectra from 501 *Apo-Ident* customers from a total of 27 007 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Myrtle oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Myrtle oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	100	0	93 774
Type C	0	2	0	91 350

The substance/substance group *Myrtle oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6630 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31836	31836	0.00	28.27
32349	31836	3.37	29.77
33879	33879	0.00	27.06

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Non-ionic hydrophil ointment
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30463-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Non-ionic hydrophil ointment; Non-ionic hydrophil ointment, conserved; Unguentum emulsificans nonionicum aquosum DAB

Special notes

When selecting the *Non-ionic hydrophil ointment* substance/substance group, the following information is displayed to the user:

- Advice for the identification of Non-ionic hydrophil ointment DAB: Improper storage of samples can result in identification problems. Please homogenises the sample and let it adapt to room temperature before the measurement.

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Non-ionic hydrophil ointment*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Non-ionic hydrop...	13138001	31001	40	20130531
Caelo	Non-ionic hydrop...	13332901	31227	40	20131104

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Non-ionic hydrop...	13333101	31451	40	1509156
Caelo	Non-ionic hydrop...	15151601	31748	40	9129
Bombastus	Non-ionic hydrop...	295368	31749	30	20150625
Caelo	Non-ionic hydrop...	161042	32088	20	20160411*
Bombastus	Non-ionic hydrop...	300481	32105	40	20160225
Bombastus	Non-ionic hydrop...	300481	32726	40	20160225
Fagron	Non-ionic hydrop...	170744	33360	40	20170413*
Caelo	Non-ionic hydrop...	172006	34010	40	20170810*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 370 spectra of 10 reference samples from the substance/substance group *Non-ionic hydrophil ointment*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 9 different batches.
- 35 174 spectra from a total of 831 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 1460 spectra from 38 reference samples of the substance/substance group *Non-ionic hydrophil ointment*.
- These include spectra of independent samples from 26 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Non-ionic hydrophil ointment	11273301	30463	40
Caelo	Non-ionic hydrophil ointment	15151601	31748 [†]	20
Bombastus	Non-ionic hydrophil ointment	295368	31749 [†]	30
Caelo	Non-ionic hydrophil ointment	153446	32077	40
Caelo	Non-ionic hydrophil ointment	153447	32078	40
Caelo	Non-ionic hydrophil ointment	153463	32084	40
Caelo	Non-ionic hydrophil ointment	161045	32085	40
Caelo	Non-ionic hydrophil ointment	161041	32087	40
Caelo	Non-ionic hydrophil ointment	161042	32088 [†]	20
Caelo	Non-ionic hydrophil ointment	153465	32106	40
Caelo	Non-ionic hydrophil ointment	160438	32107	40
Caelo	Non-ionic hydrophil ointment	161041	32108	40
Caelo	Non-ionic hydrophil ointment	161042	32109	30
Caelo	Non-ionic hydrophil ointment	161043	32110	40
Caelo	Non-ionic hydrophil ointment	161043	32111	40
Caelo	Non-ionic hydrophil ointment	161044	32112	40
Caelo	Non-ionic hydrophil ointment	161044	32113	40
Caelo	Non-ionic hydrophil ointment	161045	32114	40
Caelo	Non-ionic hydrophil ointment	161735	32417	40
Caelo	Non-ionic hydrophil ointment	161736	32418	40
Caelo	Non-ionic hydrophil ointment	161793	32419	40
Caelo	Non-ionic hydrophil ointment	161793	32420	40
Caelo	Non-ionic hydrophil ointment	161794	32421	40
Caelo	Non-ionic hydrophil ointment	162348	32488	40
Caelo	Non-ionic hydrophil ointment	162349	32489	40
Caelo	Non-ionic hydrophil ointment	162350	32490	40
Caelo	Non-ionic hydrophil ointment	170866	33361	40
Caelo	Non-ionic hydrophil ointment	170745	33362	40
Caelo	Non-ionic hydrophil ointment	170865	33363	40
Caelo	Non-ionic hydrophil ointment	170744	33364	40
Caelo	Non-ionic hydrophil ointment	170865	33365	40
Caelo	Non-ionic hydrophil ointment	170742	33366	40
Caelo	Non-ionic hydrophil ointment	170746	33367	40
Caelo	Non-ionic hydrophil ointment	170866	33368	40
Caelo	Non-ionic hydrophil ointment	170743	33369	40
Caelo	Non-ionic hydrophil ointment	171510	34009	40
Caelo	Non-ionic hydrophil ointment	172946	34011	40
Caelo	Non-ionic hydrophil ointment	181079	34154	40

- 92 414 spectra from a total of 2137 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 573 spectra from 187 *Apo-Ident* customers from 188 batches from the substance/substance group *Non-ionic hydrophil ointment*.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- These include spectra of independent samples from 180 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Non-ionic hydrophil ointment	5052015B	1
Caelo	Non-ionic hydrophil ointment	13333101	10
Caelo	Non-ionic hydrophil ointment	13333201	2
Caelo	Non-ionic hydrophil ointment	13332901	7
Caelo	Non-ionic hydrophil ointment	14068801	2
Caelo	Non-ionic hydrophil ointment	14068601	11
Caelo	Non-ionic hydrophil ointment	14259001	5
Caelo	Non-ionic hydrophil ointment	15002801	3
Caelo	Non-ionic hydrophil ointment	15002601	2
Caelo	Non-ionic hydrophil ointment	15002901	10
Caelo	Non-ionic hydrophil ointment	15007402	4
Noweda08.06.2015	Non-ionic hydrophil ointment	15027801	1
Caelo	Non-ionic hydrophil ointment	15027801	17
Caelo	Non-ionic hydrophil ointment	15105001	2
Caelo	Non-ionic hydrophil ointment	15105101	8
Sanacop	Non-ionic hydrophil ointment	15105101	1
Caelo	Non-ionic hydrophil ointment	15151801	5
Caelo	Non-ionic hydrophil ointment	15105201	1
Caelo	Non-ionic hydrophil ointment	15286601	1
Caelo	Non-ionic hydrophil ointment	15286701	4
Caelo	Non-ionic hydrophil ointment	15287001	11
Caelo	Non-ionic hydrophil ointment	15286901	4
Caelo	Non-ionic hydrophil ointment	15346401	3
Caelo	Non-ionic hydrophil ointment	15346301	1
Caelo	Non-ionic hydrophil ointment	16043801	4
Caelo	Non-ionic hydrophil ointment	16104101	11
Phoenix	Non-ionic hydrophil ointment	16104501	2
Caelo	Non-ionic hydrophil ointment	16173501	6
Caelo	Non-ionic hydrophil ointment	16179302	3
Caelo	Non-ionic hydrophil ointment	16179303	3
Caelo	Non-ionic hydrophil ointment	16179401	18
Caelo	Non-ionic hydrophil ointment	16234801	5
Bombastus	Non-ionic hydrophil ointment	17001326	1
Caelo	Non-ionic hydrophil ointment	17074201	5
Caelo	Non-ionic hydrophil ointment	17074401	2
Bombastus	Non-ionic hydrophil ointment	17086701	1
Caelo	Non-ionic hydrophil ointment	17151001	2
Caelo	Non-ionic hydrophil ointment	17294602	2
Bombastus	Non-ionic hydrophil ointment	282902	1
Bombastus	Non-ionic hydrophil ointment	282904	3
Bombastus	Non-ionic hydrophil ointment	290283	7
Caelo	Non-ionic hydrophil ointment	290283	1
Bombastus	Non-ionic hydrophil ointment	293401	2
Bombastus	Non-ionic hydrophil ointment	293402	2
Bombastus	Non-ionic hydrophil ointment	295368	5
Bombastus	Non-ionic hydrophil ointment	297633	3
Bombastus	Non-ionic hydrophil ointment	300479	4
Bombastus	Non-ionic hydrophil ointment	300481	5
Caelo	Non-ionic hydrophil ointment	300481	2
Bombastus	Non-ionic hydrophil ointment	301279	2
Bombastus	Non-ionic hydrophil ointment	303643	4
Bombastus	Non-ionic hydrophil ointment	303644	3
Bombastus	Non-ionic hydrophil ointment	308108	2
Euro OTC	Non-ionic hydrophil ointment	311098	1
Bombastus	Non-ionic hydrophil ointment	271981	3
Bombastus	Non-ionic hydrophil ointment	272063	3
Caelo	Non-ionic hydrophil ointment	11357601	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Non-ionic hydrophil ointment	11357502	1
Caelo	Non-ionic hydrophil ointment	12086901	3
Caelo	Non-ionic hydrophil ointment	12001602	1
Caelo	Non-ionic hydrophil ointment	12087001	4
Caelo	Non-ionic hydrophil ointment	12155801	2
Caelo	Non-ionic hydrophil ointment	12156101	2
phönix 22.01.2013	Non-ionic hydrophil ointment	12156101	1
Fagron	Non-ionic hydrophil ointment	12223301	3
Caelo	Non-ionic hydrophil ointment	12223501	3
Caelo	Non-ionic hydrophil ointment	12353001	1
Caelo	Non-ionic hydrophil ointment	13045401	4
Caelo	Non-ionic hydrophil ointment	12352801	7
Fagron	Non-ionic hydrophil ointment	13045301	1
Caelo	Non-ionic hydrophil ointment	13045501	3
Caelo	Non-ionic hydrophil ointment	13138001	16
Caelo	Non-ionic hydrophil ointment	13138101	1
Caelo	Non-ionic hydrophil ointment	13263301	2
Caelo	Non-ionic hydrophil ointment	13138201	4
Caelo	Non-ionic hydrophil ointment	13138101	8
Caelo	Non-ionic hydrophil ointment	13263501	7
Caelo	Non-ionic hydrophil ointment	13333401	5
Noweda	Non-ionic hydrophil ointment	1812U-02543	1
Caelo	Non-ionic hydrophil ointment	21111203	1
Bombastus	Non-ionic hydrophil ointment	280909	3
Bombastus	Non-ionic hydrophil ointment	6861335	1
GEHE 10.03.2014	Non-ionic hydrophil ointment	13137901	1
Caelo	Non-ionic hydrophil ointment	13263801	2
Caelo	Non-ionic hydrophil ointment	1333101	1
Berkel AHK	Non-ionic hydrophil ointment	13332901	3
Caelo	Non-ionic hydrophil ointment	1326501	1
VDL	Non-ionic hydrophil ointment	13332901	1
Caelo	Non-ionic hydrophil ointment	13333001	3
Berkel AHK	Non-ionic hydrophil ointment	13333101	2
Caelo	Non-ionic hydrophil ointment	13333301	3
Caelo	Non-ionic hydrophil ointment	13333501	6
Caelo	Non-ionic hydrophil ointment	14068401	1
Sanacorp	Non-ionic hydrophil ointment	14068601	1
Caelo	Non-ionic hydrophil ointment	14068701	2
Berkel AHK	Non-ionic hydrophil ointment	14068601	2
Caelo	Non-ionic hydrophil ointment	14068901	2
Caelo	Non-ionic hydrophil ointment	14069001	2
AHD	Non-ionic hydrophil ointment	14069001	3
Caelo	Non-ionic hydrophil ointment	14221301	1
Caelo	Non-ionic hydrophil ointment	14221401	6
GEHE 23.09.2014	Non-ionic hydrophil ointment	14221301	1
Caelo	Non-ionic hydrophil ointment	14221501	2
Caelo	Non-ionic hydrophil ointment	14258901	6
Caelo	Non-ionic hydrophil ointment	14259101	3
Caelo	Non-ionic hydrophil ointment	14259201	2
Sanacorp	Non-ionic hydrophil ointment	15002901	1
Caelo	Non-ionic hydrophil ointment	15027601	1
Caelo	Non-ionic hydrophil ointment	15105301	1
Caelo	Non-ionic hydrophil ointment	15151601	5
Caelo	Non-ionic hydrophil ointment	15286801	2
Caelo	Non-ionic hydrophil ointment	15344701	1
Caesar/Noweda	Non-ionic hydrophil ointment	15346401	1
Bombastus	Non-ionic hydrophil ointment	160905UEN	1
Caelo	Non-ionic hydrophil ointment	16104201	1
Gehe	Non-ionic hydrophil ointment	16104401	1
Bombastus	Non-ionic hydrophil ointment	161220NHC/301280	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Non-ionic hydrophil ointment	16173401	3
Noweda	Non-ionic hydrophil ointment	16179302	1
Caelo	Non-ionic hydrophil ointment	16179301	1
Caelo	Non-ionic hydrophil ointment	16179402	1
Caelo	Non-ionic hydrophil ointment	16179603	1
Noweda	Non-ionic hydrophil ointment	16234801	1
Caelo	Non-ionic hydrophil ointment	16234901	2
Caelo	Non-ionic hydrophil ointment	16235001	1
Caelo	Non-ionic hydrophil ointment	170328UENA303640	1
Caelo	Non-ionic hydrophil ointment	17074301	4
Caelo	Non-ionic hydrophil ointment	17074501	1
Caelo	Non-ionic hydrophil ointment	17086501	3
Caelo	Non-ionic hydrophil ointment	17086701	1
Bombastus	Non-ionic hydrophil ointment	170905NHC/308112	1
Bombastus	Non-ionic hydrophil ointment	171128NHC/308111	1
Bombastus	Non-ionic hydrophil ointment	171114NHC/308111	1
Caelo	Non-ionic hydrophil ointment	17112014A	1
Caelo	Non-ionic hydrophil ointment	17151101	1
Noweda	Non-ionic hydrophil ointment	17994601	1
Bombastus	Non-ionic hydrophil ointment		1
Bombastus	Non-ionic hydrophil ointment	180308NHC/308108	1
Caelo	Non-ionic hydrophil ointment	19081505	1
Caelo	Non-ionic hydrophil ointment	20061406	1
Bombastus	Non-ionic hydrophil ointment	286002	2
Bombastus	Non-ionic hydrophil ointment	286003	5
Bombastus	Non-ionic hydrophil ointment	286001	1
bom	Non-ionic hydrophil ointment	286003	1
Bombastus	Non-ionic hydrophil ointment	288042	3
Bombastus	Non-ionic hydrophil ointment	288043	6
Sanacorp	Non-ionic hydrophil ointment	288043	1
Bombastus	Non-ionic hydrophil ointment	288044	6
Caelo	Non-ionic hydrophil ointment	288044	1
Bombastus	Non-ionic hydrophil ointment	290281	1
Fagron	Non-ionic hydrophil ointment	290282	1
Bombastus	Non-ionic hydrophil ointment	293403	2
Bombastus	Non-ionic hydrophil ointment	295367	1
Bombastus	Non-ionic hydrophil ointment	296801	4
Bombastus	Non-ionic hydrophil ointment	296803	2
Bombastus	Non-ionic hydrophil ointment	296802	1
Caelo	Non-ionic hydrophil ointment	297634	1
Bombastus	Non-ionic hydrophil ointment	297634	1
Bombastus	Non-ionic hydrophil ointment	300480	1
Bombastus	Non-ionic hydrophil ointment	301280	4
Bombastus	Non-ionic hydrophil ointment	303640	5
Bombastus	Non-ionic hydrophil ointment	303642	1
Bombastus	Non-ionic hydrophil ointment	306522	4
Bombastus	Non-ionic hydrophil ointment	308109	1
Bombastus	Non-ionic hydrophil ointment	308111	1
Bombastus	Non-ionic hydrophil ointment	308112	1
Bombastus	Non-ionic hydrophil ointment	258174	1
Bombastus	Non-ionic hydrophil ointment	266631	2
Euro OTC	Non-ionic hydrophil ointment	266631	1
Bombastus	Non-ionic hydrophil ointment	267625	3
Bombastus	Non-ionic hydrophil ointment	268754	2
Bombastus	Non-ionic hydrophil ointment	271976	1
Sanacorp	Non-ionic hydrophil ointment	271980	1
Bombastus	Non-ionic hydrophil ointment	271977	1
Noweda	Non-ionic hydrophil ointment	271980	2
Sanacorp	Non-ionic hydrophil ointment	271981	1
Bombastus, Noweda	Non-ionic hydrophil ointment	272063	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Non-ionic hydrophil ointment	2122013c	1
Caelo	Non-ionic hydrophil ointment	513Q-02543	1
Caelo	Non-ionic hydrophil ointment	6061305	1
Caelo	Non-ionic hydrophil ointment	11273001	1
Caelo	Non-ionic hydrophil ointment	11273303	2
Caelo	Non-ionic hydrophil ointment	11273401	1
Caelo	Non-ionic hydrophil ointment	11357301	1
Caelo	Non-ionic hydrophil ointment	11357801	1
Caelo	Non-ionic hydrophil ointment	12086802	1
Phönix 25.09.2012	Non-ionic hydrophil ointment	12087001	1
Caelo	Non-ionic hydrophil ointment	12086902	1
Caelo	Non-ionic hydrophil ointment	12155901	5
Caelo	Non-ionic hydrophil ointment	12223301	1
Caelo	Non-ionic hydrophil ointment	12223602	1
Caelo	Non-ionic hydrophil ointment	12223601	3
Caelo	Non-ionic hydrophil ointment	12352901	3
Caesar &Loretz GmbH	Non-ionic hydrophil ointment	12352901	1
Caelo	Non-ionic hydrophil ointment	130450401	1
Caelo	Non-ionic hydrophil ointment	13045201	3
Sanacorp	Non-ionic hydrophil ointment	13045401	1
Caelo	Non-ionic hydrophil ointment	13045601	3
Hustadt Apotheke	Non-ionic hydrophil ointment	13045601	1
Caelo	Non-ionic hydrophil ointment	13137801	1
Caelo	Non-ionic hydrophil ointment	13137901	1
Martinus-Apotheke	Non-ionic hydrophil ointment	13138001	1
Noweda	Non-ionic hydrophil ointment	13138101	1
Post-Apotheke	Non-ionic hydrophil ointment	13138301	1
Caelo	Non-ionic hydrophil ointment	1318001	1
Caelo	Non-ionic hydrophil ointment	13263401	1
Caelo	Non-ionic hydrophil ointment	13263601	3
Caelo	Non-ionic hydrophil ointment	13263701	4
Post-Apotheke	Non-ionic hydrophil ointment	13263701	2
Caelo	Non-ionic hydrophil ointment	13263902	3
Caelo	Non-ionic hydrophil ointment	21081202	1
Bombastus	Non-ionic hydrophil ointment	271973	1
Bombastus	Non-ionic hydrophil ointment	275256	1
Bombastus	Non-ionic hydrophil ointment	275259	2
Bombastus	Non-ionic hydrophil ointment	275260	2
Bombastus	Non-ionic hydrophil ointment	278694	2
Caelo	Non-ionic hydrophil ointment	278694	1
Bombastus	Non-ionic hydrophil ointment	282903	2
Caelo	Non-ionic hydrophil ointment		1
Caelo	Non-ionic hydrophil ointment	51121-02543	1
Caelo	Non-ionic hydrophil ointment	4812I-02543	1
Caelo	Non-ionic hydrophil ointment	CHB13263501	1
Caelo	Non-ionic hydrophil ointment	11273302	1
Noweda	Non-ionic hydrophil ointment	11357301	1
Caelo	Non-ionic hydrophil ointment	11357501	1
Sanacorp	Non-ionic hydrophil ointment	13138001	1

- 90 779 spectra from 501 *Apo-Ident* customers from a total of 26 826 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Non-ionic hydrophil ointment* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Non-ionic hydrophil ointment* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	370	0	35 174
Type B	0	1406	54	92 414
Type C	0	485	88	90 779

The substance/substance group *Non-ionic hydrophil ointment* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.3784 %)
Type B	100.0000 % (> 99.9761 %)	96.3014 % (> 96.0959 %)
Type C	100.0000 % (> 99.6586 %)	84.6422 % (> 84.1187 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31001	31748	9.60	17.08
31227	31748	7.51	16.53
31451	31451	0.00	17.34
31748	31748	0.00	13.54
31749	31748	4.25	16.79
32088	32088	0.00	13.01
32105	31748	5.27	16.12
32726	33360	7.52	19.20
33360	33360	0.00	19.17
34010	34010	0.00	11.55

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Octyldodecanol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30867-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Octyldodecanol; Octyldodecanolum

Special notes

When selecting the *Octyldodecanol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Octyldodecanol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Octyldodecanol	12165509	30867	40	20120711
Caelo	Octyldodecanol	13360010	31330	40	AR-15-FG-000314-01
Caelo	Octyldodecanol	15135212	31753	40	20150506
Caelo	Octyldodecanol	16095507	32958	40	20160421

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Octyldodecanol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 176 spectra from 6 reference samples of the substance/substance group *Octyldodecanol*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Octyldodecanol	15135212	31753 [†]	20
Caelo	Octyldodecanol	163127	33093	40
Caelo	Octyldodecanol	171381	34052	40
Caelo	Octyldodecanol	181120	34220	40
Caelo	Octyldodecanol	16095506	not required	24
Fagron	Octyldodecanol	16G20-B99-324356	not required	12

- 93 698 spectra from a total of 2161 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1204 spectra from 307 *Apo-Ident* customers from 174 batches from the substance/substance group *Octyldodecanol*.
- These include spectra of independent samples from 168 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Octyldodecanol	14371706	17
Caelo	Octyldodecanol	14371711	25
Caelo	Octyldodecanol	13222011	8
Phoenix	Octyldodecanol	13360006	1
Caelo	Octyldodecanol	13360008	14
Caelo	Octyldodecanol	13360010	15
Caelo	Octyldodecanol	14130205	14
Caelo	Octyldodecanol	14130206	16
Caelo	Octyldodecanol	13360011	14
Caelo	Octyldodecanol	14130207	30
Caelo	Octyldodecanol	14158410	9
Caelo	Octyldodecanol	2613M-02875	2
Caelo	Octyldodecanol	10281907	1
Caelo	Octyldodecanol	3314E-02875	1
Gehe	Octyldodecanol	11060308	1
Caelo	Octyldodecanol	11060312	1
Caelo	Octyldodecanol	11060319	3
Caelo	Octyldodecanol	1136017	1
Caelo	Octyldodecanol	11369003	8
Caelo	Octyldodecanol	11369004	4
Caelo	Octyldodecanol	11369012	2
Gehe	Octyldodecanol	11369012	1
Caelo	Octyldodecanol	11369014	2
Caelo	Octyldodecanol	11369016	8
Caelo	Octyldodecanol	11369017	6
phönix, 20.11.2012	Octyldodecanol	11369017	2
Caelo	Octyldodecanol	11369021	6
Gehe	Octyldodecanol	11369021	1
Gehe	Octyldodecanol	12165506	1
Caelo	Octyldodecanol	12165506	4
Noweda	Octyldodecanol	11369021	1
Caelo	Octyldodecanol	12091304	1
Caelo	Octyldodecanol	12165509	6
Anzag	Octyldodecanol	12165509	1
Caelo	Octyldodecanol	12165511	11
phönix, 05.12.2012	Octyldodecanol	12165511	1
Caelo	Octyldodecanol	12165514	4
Sanacorp	Octyldodecanol	12165514	2
Caelo	Octyldodecanol	12165515	10
Noweda	Octyldodecanol	12165515	1
Anzag, 20.03.2013	Octyldodecanol	12165515	2
Gehe	Octyldodecanol	12165515	2
Caelo	Octyldodecanol	12165517	2
Caelo	Octyldodecanol	12165519	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Octyldodecanol	12165519	1
Caelo	Octyldodecanol	12165520	12
Caelo	Octyldodecanol	12165521	12
Phönix	Octyldodecanol	12165521	3
Caelo	Octyldodecanol	13044002	10
Phönix	Octyldodecanol	13044005	1
Beiersdorf/Anzag	Octyldodecanol	13044007	1
Caelo	Octyldodecanol	13044010	12
Caelo	Octyldodecanol	13044013	20
Alliance healthcare	Octyldodecanol	13044013	1
Caelo	Octyldodecanol	13044016	21
Caelo	Octyldodecanol	13222003	20
Caelo	Octyldodecanol	13044015	8
Phönix	Octyldodecanol	13222008	1
Caelo	Octyldodecanol	13222008	18
Caelo	Octyldodecanol	13360012	14
Caelo	Octyldodecanol	15135206	31
Caelo	Octyldodecanol	14371714	4
Caelo	Octyldodecanol	15135203	14
Caelo	Octyldodecanol	13360002	14
Caelo	Octyldodecanol	15135211	9
Caelo	Octyldodecanol	15135212	10
Caelo	Octyldodecanol	15383004	15
Caelo	Octyldodecanol	15135213	59
Caelo	Octyldodecanol	15383005	39
Caelo	Octyldodecanol	16095502	10
Ichthyol/Phönix	Octyldodecanol	16095503	1
Caelo	Octyldodecanol	26071208	1
Fagron	Octyldodecanol	200714112-2	1
Caelo	Octyldodecanol	82695089	1
Caelo	Octyldodecanol		2
Caelo	Octyldodecanol	B11060317	1
Fagron	Octyldodecanol	10598578	1
Caelo	Octyldodecanol	516A-02875	1
Caelo	Octyldodecanol	11052015A	1
Caelo	Octyldodecanol	8012015B	1
GEHE	Octyldodecanol	1098685	1
Nachprüfung	Octyldodecanol	13044015	2
Caelo	Octyldodecanol	1314I-02875	1
Caelo	Octyldodecanol	13222002	11
Caelo	Octyldodecanol	13222013	1
Caelo	Octyldodecanol	133360010	1
Caesar/Noweda	Octyldodecanol	13222010	1
Caelo	Octyldodecanol	13360001	3
Noweda	Octyldodecanol	13360003	1
Caelo	Octyldodecanol	13360003	8
Caelo	Octyldodecanol	13360005	6
Klenk	Octyldodecanol	13360003	1
Caesar & Loretz	Octyldodecanol	13360003	1
Caelo	Octyldodecanol	13360006	12
VDL 30.04.14_5.77 euro	Octyldodecanol	13360006	1
Gehe	Octyldodecanol	13360006	1
Caelo	Octyldodecanol	13360009	2
Caesar & Loretz GmbH	Octyldodecanol	13360010	1
AHD 30.5 preis:6,07	Octyldodecanol	13360011//05R06	1
PHOENIX 30.06.2014	Octyldodecanol	13360011	1
Noweda	Octyldodecanol	13360012	1
Fagron	Octyldodecanol	13L02-B05	2
Sanacorp	Octyldodecanol	13102-b05-288974	1
Caelo	Octyldodecanol	14130202	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Octyldodecanol	14130204	3
Fagron	Octyldodecanol	13L02-B05-288975	1
Caeo	Octyldodecanol	14130205	1
Phönix	Octyldodecanol	14130206	1
sanacorp	Octyldodecanol	14130207	2
Phoenix	Octyldodecanol	14130207	1
Caesar & Loretz	Octyldodecanol	14130207	1
Caelo	Octyldodecanol	14135206	1
Caelo	Octyldodecanol	14158401	1
Caelo	Octyldodecanol	14158402	15
VDL; 5,77EUR;21.01.2015	Octyldodecanol	14158402	1
Caelo	Octyldodecanol	14158407	1
Caelo	Octyldodecanol	14158411	2
Phoenix	Octyldodecanol	14158411	1
Caelo	Octyldodecanol	1431707	1
Caelo	Octyldodecanol	143710706	1
Bombastus	Octyldodecanol	14371703	1
Caelo	Octyldodecanol	14371703	1
Caelo	Octyldodecanol	14371705	5
sanacorp	Octyldodecanol	14371706	1
Caelo	Octyldodecanol	14371707	4
Caelo	Octyldodecanol	14371710	2
VDL,5,77EUR; 27.04.15	Octyldodecanol	14371711	2
Sanacorp	Octyldodecanol	14371711	1
Anzag/Henry Lamotte	Octyldodecanol	14371711	1
Noweda 10.07.2015	Octyldodecanol	14371714	1
Caelo	Octyldodecanol	14371716	9
Sanacorp	Octyldodecanol	14371716	1
Fagron	Octyldodecanol	14E06-B02	4
Fagron	Octyldodecanol	14e06-b02-293745	1
Fagron	Octyldodecanol	14E06-B02-293746	1
Fagron	Octyldodecanol	14E06-B02-30170	1
Fagron	Octyldodecanol	14E06-B02-293745	1
Fagron	Octyldodecanol	14E06-B02-301705	2
Fagron	Octyldodecanol	14E06-B02-301706	1
Caelo	Octyldodecanol	1513513	1
Caelo	Octyldodecanol	15092014B	1
Caelo	Octyldodecanol	15135205	1
caesar&loretz,Phönix	Octyldodecanol	15135206	1
Caesar & Loretz GmbH	Octyldodecanol	15135206	1
Caelo	Octyldodecanol	15135207	7
Caesar/Noweda	Octyldodecanol	15135207	1
Noweda	Octyldodecanol	15135213	1
Fagron	Octyldodecanol	15135213	1
Caelo	Octyldodecanol	1538005	2
Caelo	Octyldodecanol	1538004	1
Caelo	Octyldodecanol	1538003	1
Caelo	Octyldodecanol	153830004	1
Phönix	Octyldodecanol	1538009	1
Caelo	Octyldodecanol	1538009	1
Caelo	Octyldodecanol	15383002	1
Caelo	Octyldodecanol	15383003	6
Fagron	Octyldodecanol	15383004	1
Euro OTC	Octyldodecanol	15383004	1
Caesar&Loretz GmbH	Octyldodecanol	15383004	1
VDL;15.3.16;5,75EUR	Octyldodecanol	15383005	1
Hustadt-Apotheke	Octyldodecanol	15383005	1
Caolo	Octyldodecanol	15383005	1
Caeol/Phoenix	Octyldodecanol	15383005	1
Caelo	Octyldodecanol	15383008	21

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Octyldodecanol	15383008	1
Caeolo/Noweda	Octyldodecanol	15383008	1
Caelo	Octyldodecanol	15383009	19
Noweda	Octyldodecanol	15383009	1
Phoenix	Octyldodecanol	15383009	1
Caelo	Octyldodecanol	1583009	1
Caelo	Octyldodecanol	1583004	1
Caelo	Octyldodecanol	16095501	1
Fagron	Octyldodecanol	15E06-B02-30170	1
Caesar&Loretz	Octyldodecanol	16095502	1
Caelo	Octyldodecanol	16095503	20
Caesar & Loretz/AHC	Octyldodecanol	16095503	1
Phoenix	Octyldodecanol	16095503	1
Caesar & Loretz/Noweda	Octyldodecanol	16095503	2
Gehe	Octyldodecanol	16095503	1
Caelo	Octyldodecanol	16095506	8
Phönix	Octyldodecanol	16095506	1
Caelo	Octyldodecanol	16095507	34
Bombastus	Octyldodecanol	16095506	1
GEHE	Octyldodecanol	16095507	1
Anzag	Octyldodecanol	16095507	1
Caesar&Loretz GmbH,09.0...	Octyldodecanol	16095507	1
Caelo	Octyldodecanol	16095508	1
Caelo	Octyldodecanol	16095511	21
Phoenix WE 22.03.17 EK...	Octyldodecanol	16095511	1
Caelo	Octyldodecanol	16312701	26
Noweda	Octyldodecanol	16095511	1
Fagron	Octyldodecanol	16312701	1
Euro OTC	Octyldodecanol	16312701	1
Gehe	Octyldodecanol	16312701	2
Phoenix	Octyldodecanol	16312701	1
Krieger	Octyldodecanol	16312701	1
AHD	Octyldodecanol	16312701	1
GEHE	Octyldodecanol	16312701	1
Caelo	Octyldodecanol	16312704	1
Caelo	Octyldodecanol	16312706	12
Caelo	Octyldodecanol	16312707	7
Caelo	Octyldodecanol	16312710	4
Caeolo / Phönix	Octyldodecanol	16312707	1
Caelo	Octyldodecanol	17077701	2
Caelo	Octyldodecanol	17082015C	1
Caelo	Octyldodecanol	17138102	12
Fagron	Octyldodecanol	17138102	1
Cealo/Noweda	Octyldodecanol	17138102	1
Sanacorp	Octyldodecanol	17138102	1
Gehe	Octyldodecanol	17138102	1
Caelo	Octyldodecanol	17138103	6
Sanacorp WE:09.11.17 10...	Octyldodecanol	17138105	1
Noweda	Octyldodecanol	17138105	1
Caelo	Octyldodecanol	17138105	4
gehe	Octyldodecanol	17138109	1
Sanacorp WE:21.12.17 E...	Octyldodecanol	17138109	1
Caelo	Octyldodecanol	17138109	1
Caelo	Octyldodecanol	17273001	5
Phönix W.E:27.02.18 E...	Octyldodecanol	17273001	1
Caelo	Octyldodecanol	17273002	11
GEHE	Octyldodecanol	17273002	1
Caelo	Octyldodecanol	1738102	1
Caelo	Octyldodecanol	1801190D/17138109	1
Caelo	Octyldodecanol	1802080D/17138109	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Octyldodecanol	17273006	1
Caelo	Octyldodecanol	18112001	1
Caelo	Octyldodecanol	24062014B	1
Caelo	Octyldodecanol	31101416	1
Caelo	Octyldodecanol	4214U-02875	1
Caelo	Octyldodecanol	11060305	1
Caeol	Octyldodecanol	6022014B	1
Gehe	Octyldodecanol	1136017	1
Caelo	Octyldodecanol	11012014A	1
Spangropharm	Octyldodecanol	11060317	1
Phönix 08.11.12	Octyldodecanol	11369021	1
Fargon	Octyldodecanol	11C14-N14	1
Fargon	Octyldodecanol	11J13-N05	5
Caelo	Octyldodecanol	12044009	1
Caelo	Octyldodecanol	12111303	1
Caelo	Octyldodecanol	12011401	1
Sanacorp	Octyldodecanol	12142726	1
Noweda	Octyldodecanol	1216551	1
Kehr	Octyldodecanol	12165519	1
Anzag	Octyldodecanol	12165519	1
AHD 09.04.2013	Octyldodecanol	12165521	1
Caelo	Octyldodecanol	1265515	2
Fargon	Octyldodecanol	12C21-N03	3
Fargon	Octyldodecanol	12J04-N01	3
Fargon	Octyldodecanol	12C31-N03	1
Phoenix/Fargon	Octyldodecanol	12J04-N01	1
Phönix Gotha	Octyldodecanol	13044002	1
Caelo	Octyldodecanol	13044004	1
Noweda	Octyldodecanol	13044005	1
Caelo	Octyldodecanol	13044005	14
Sanacorp	Octyldodecanol	13044005	2
Caelo	Octyldodecanol	13044007	8
Caesar/Lorenz	Octyldodecanol	13044009	1
Caelo	Octyldodecanol	13044009	14
Gehe 27.06.2013	Octyldodecanol	13044009	1
Caesar & Loretz/ Noweda	Octyldodecanol	13044009	1
Caesar/Noweda	Octyldodecanol	13044009	1
Noweda	Octyldodecanol	13044013	2
Gehe	Octyldodecanol	13044013	1
Euro OTC	Octyldodecanol	13044015	1
Noweda	Octyldodecanol	13044016	1
AHD, 09.10.2013	Octyldodecanol	13044016	1
Caelo	Octyldodecanol	13092013A	1
Caesar&lorenz,Phoenix	Octyldodecanol	13044016	1
Caelo	Octyldodecanol	1322005	1
Sanacorp	Octyldodecanol	13222002	1
Noweda	Octyldodecanol	13222003	1
Sanacorp	Octyldodecanol	13222003	4
Caelo	Octyldodecanol	13222005	4
Noweda	Octyldodecanol	13222008	1
Sanacorp	Octyldodecanol	13222008	1
Caesar & Loretz GmbH	Octyldodecanol	13222008	1
Cealo / Noweda	Octyldodecanol	13222008	1
Caelo	Octyldodecanol	13222010	1
Noweda	Octyldodecanol	13360002	1
Sanacorp	Octyldodecanol	13360006	1
Caesar&lorenz,Phoenix	Octyldodecanol	13360006	1
Krieger	Octyldodecanol	13360006	1
Fargon	Octyldodecanol	13E07-N03	3
Sanacorp	Octyldodecanol	13E07-N07	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar&lorenz,Phoenix	Octyldodecanol	14371711	1
Caelo	Octyldodecanol	2513A-02875	2
Caelo	Octyldodecanol	2xa13360002	1
Caelo	Octyldodecanol	4113A2875	1
Sanacorp	Octyldodecanol	AKD31	1
Fagron	Octyldodecanol	CE11350022	2
Caelo	Octyldodecanol	CE11350022	2
Caelo	Octyldodecanol	11060317	1
Caelo	Octyldodecanol	1098722	1
Sanacorp	Octyldodecanol	13044002	1
A.H.D.	Octyldodecanol	13044007	1
Caelo	Octyldodecanol	13044040	1
Caelo	Octyldodecanol	13222001	1
Gehe	Octyldodecanol	13360003	1
Gehe	Octyldodecanol	13360012	1
Noweda	Octyldodecanol	21022171/7	1

- 90 148 spectra from 501 *Apo-Ident* customers from a total of 26 840 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Octyldodecanol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Octyldodecanol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	176	0	93 698
Type C	0	1192	12	90 148

The substance/substance group *Octyldodecanol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.5909 %)
Type C	100.0000 % (> 99.6586 %)	99.0033 % (> 98.7542 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30867	31330	3.94	54.20
31330	31330	0.00	52.77
31753	31330	12.40	54.89
32958	31330	7.93	50.90

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Oleic acid**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 32471-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Oleic acid; Acidum oleicum; Acidum oleinicum

Special notes

When selecting the *Oleic acid* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Oleic acid*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Oleic acid	15298701	32471	40	1702061
Caelo	Oleic acid	163091	33076	40	20161212*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Oleic acid*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Oleic acid*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Oleic acid	15298701	32475	40
Caelo	Oleic acid	15298701	32747	40
Caelo	Oleic acid	172741	33898	40

- 93 754 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 0 spectra from 0 *Apo-Ident* customers from 0 batches from the substance/substance group *Oleic acid*.
- These include spectra of independent samples from 0 batches, of which no spectra have been used for the generation of the database.
- 91 352 spectra from 501 *Apo-Ident* customers from a total of 27 009 batches from a further 334 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Oleic acid* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Oleic acid* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	120	0	93 754
Type C	0	0	0	91 352

The substance/substance group *Oleic acid* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6586 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
32471	32471	0.00	56.82
33076	33076	0.00	53.79

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Oregano oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30669-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Oregano oil; Oleum origani cretici

Special notes

When selecting the *Oregano oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Oregano oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Oregano oil	11247605	30669	40	20110829
Caelo	Oregano oil	13354002	31196	40	1502253
Caelo	Oregano oil	13354007	31302	40	20131028
Caelo	Oregano oil	13354006	31406	40	20131028
Caelo	Oregano oil	14242905	31788	40	20140904

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Oregano oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 344 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Oregano oil*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Oregano oil	14242905	31788 [†]	20
Caelo	Oregano oil	16011001	32070	60
Caelo	Oregano oil	16011001	33469	40

- 93 754 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 4 spectra from 4 *Apo-Ident* customers from 2 batches from the substance/substance group *Oregano oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Noweda/Taoasis	Oregano oil	43878A-123313	1
Taoasis	Oregano oil	43878A-123313	1
Caelo	Oregano oil	12347604	2

- 91 348 spectra from 501 *Apo-Ident* customers from a total of 27 007 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Oregano oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Oregano oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	120	0	93 754
Type C	0	4	0	91 348

The substance/substance group *Oregano oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6608 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30669	31196	6.56	56.10
31196	31196	0.00	50.22
31302	31196	2.94	51.38
31406	31196	6.52	52.76
31788	31196	7.93	49.85

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Paraffin classification
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30874-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Paraffin classification; Emulsifying hydrophobic base gel DAC; Eye ointment, emulsifying DAC; Eye paraffin; Gelatum basalis hydrophobicum; Gelatum basalis hydrophobicum emulsificans; Hydrophilic ointment DAB; Hydrophobic base gel DAC; Lanae alcoholum unguentum; Oculentum emulsificans DAC; Oculentum simplex DAC; Paraffin, liquid; Paraffin, white; Paraffin, yellow; Paraffinum; Unguentum alcoholum lanae; Unguentum alcoholum lanae SR; Unguentum Cordes®; Unguentum emulsificans DAB; Unguentum ophthalmicum emulsificans DAC; Vaseline, white; Vaseline, yellow; Vaselinum album; Vaselinum album pro ophthalmicum; Vaselinum flavum; Wool alcohol ointment DAB; Wool alcohol ointment SR with white paraffin; Wool alcohol ointment SR/DAC

Special notes

When selecting the *Paraffin classification* substance/substance group, the following information is displayed to the user:

- Advice for the identification of paraffin: Paraffin, liquid and paraffin, light liquid can not be conclusively distinguished by spectral analysis. An additional analysis of identity may be required.

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Paraffin classification*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Paraffin, liquid	12102101	30874	40	20120416
Caelo	Paraffin, liquid	12324307	30877	40	20121031
Caelo	Hydrophobic base...	12229003	30932	40	1405121
Caelo	Emulsifying hydr...	12343501	30933	40	1405098
Caelo	Paraffin, white	12374502	31009	40	20130102
Caelo	Paraffin, white	12032413	31010	40	20120203
Caelo	Hydrophobic base...	13299604	31168	40	20130909
apomix	Eye ointment, em...	15B241013	31174	39	1405331
Caelo	Wool alcohol oin...	12244302	31228	70	AR-15-FG-001376-01
Caelo	Wool alcohol oin...	14032501	31272	40	1507263
Hansen & Rose...	Hydrophobic base...	1023472	31292	40	20140328
Caelo	Paraffin, liquid	13248007	31351	40	AR-15-FG-002534-01
Caelo	Paraffin, liquid	13390602	31358	40	AR-15-FG-002530-01
Caelo	Paraffin, white	13341601	31437	40	20131108
Caelo	Wool alcohol oin...	14067001	31441	40	20140307
Caelo	Eye paraffin	13338710	31445	40	20131108
Caelo	Hydrophobic base...	14044601	31446	40	AR-15-FG-007940-01
Caelo	Hydrophilic oint...	14024802	31456	40	AR-14-FG-012697-01
Caelo	Paraffin, white	13268421	31468	40	20130823
Caelo	Paraffin, yellow	12321908	31473	40	AR-15-FG-002511-01
Euro OTC	Wool alcohol oin...	1211017-01	31522	40	AR-14-FG-012695-01
Euro OTC	Paraffin, yellow	1401017-01	31541	40	20140210
Caelo	Hydrophilic oint...	14024802	31567	40	20140129
Caelo	Emulsifying hydr...	14276401	31570	40	AR-15-FG-014916-01
Euro OTC	Paraffin, liquid	1409007-01	31587	40	20140929
Euro OTC	Paraffin, white	1405014-01	31590	40	20140702
Bombastus	Paraffin, white	290271	31601	40	20140929
Ichthyol	Unguentum Cordes®	13A067	31602	40	1507406
Caelo	Eye paraffin	13338711	31610	40	AR-15-FG-008217-01
Caelo	Hydrophilic oint...	13245302	31611	40	20131008
Fagron	Hydrophilic oint...	14G08-T08-013954	31612	40	20140829
Euro OTC	Paraffin, liquid	1506021	31702	40	0000010380
Euro OTC	Paraffin, white	1509019	31736	60	20150929*
Euro OTC	Paraffin, liquid	1509011	31740	40	20150918*
apomix	Wool alcohol oin...	02A070515	31744	40	AR-16-FG-004495-01
Ichthyol	Unguentum Cordes®	15A137	31807	40	20150804
Euro OTC	Paraffin, liquid	1511030	31812	40	20151208*
Euro OTC	Paraffin, liquid	1511029	31816	40	20151208*
Caelo	Hydrophilic oint...	154164	31909	40	20151215*
Caelo	Hydrophilic oint...	154165	31910	40	20151215*
Caelo	Paraffin, white	153276	31916	40	20151001*
Caelo	Paraffin, white	153368	31917	40	20151012*
Caelo	Emulsifying hydr...	153311	31920	40	20151006*
Caelo	Wool alcohol oin...	153347	31921	40	20150930*
Caelo	Paraffin, yellow	153304	31940	40	20151001*
Caelo	Paraffin, yellow	153305	31941	40	20151001*
apomix	Wool alcohol oin...	01A180215	31949	40	20150223
Ichthyol	Unguentum Cordes®	15A150	31954	40	20150930
Caelo	Wool alcohol oin...	153421	31999	40	20151012*
Caelo	Wool alcohol oin...	153562	32017	40	20151028*
Euro OTC	Paraffin, liquid	1604009	32279	40	20160511*
Caelo	Paraffin, white	162274	32495	40	20160804*
Caelo	Wool alcohol oin...	162513	32517	40	20160907*
Euro OTC	Paraffin, liquid	1607039	32704	40	20160826*
Caelo	Eye paraffin	162718	32900	40	20161020*
apomix	Eye ointment, em...	06A090516	32946	40	20160517
apomix	Eye ointment, em...	10A041016	32947	40	20161017
Caelo	Emulsifying hydr...	170323	33133	40	20170206*
Fagron	Paraffin, liquid	16E30-B05-327577	33274	40	20160614
Caelo	Emulsifying hydr...	171137	33396	40	20170426*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Paraffin, liquid	15B18-B01-333364	33729	40	20150302

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 2489 spectra of 61 reference samples from the substance/substance group *Paraffin classification*. These samples are listed above in the section *calibration samples*. The reference samples originate from 60 different batches.
- 33 055 spectra from a total of 780 batches from further 200 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 5609 spectra from 149 reference samples of the substance/substance group *Paraffin classification*.
- These include spectra of independent samples from 110 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Paraffin, white	90292119	30006	39
Caelo	Paraffin, white	10094704	30164	40
Ichthyol	Unguentum Cordes®	07A050	30218	50
Euro OTC	Wool alcohol ointment DAB	1007002-01	30269	40
Caelo	Emulsifying hydrophobic base...	10307401	30318	40
Ichthyol	Unguentum Cordes®	11A150	30375	30
Caelo	Hydrophilic ointment DAB	11211502	30393	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Ichthyol	Unguentum Cordes [®]	11A070	30403	40
Ichthyol	Unguentum Cordes [®]	11A161	30424	40
Caelo	Eye paraffin	09274706	30434	40
Caelo	Paraffin, white	11374906	30466	70
Ichthyol	Unguentum Cordes [®]	11A139	30477	40
apomix	Eye ointment, emulsifying DAC	03B270212	30493	40
Caelo	Wool alcohol ointment SR wit...	11278101	30497	40
Ichthyol	Unguentum Cordes [®]	11A199	30587	20
Caelo	Hydrophilic ointment DAB	11211601	30609	40
Ichthyol	Unguentum Cordes [®]	11A202	30626	40
Caelo	Paraffin, liquid	13248007	31351 [†]	20
Caelo	Paraffin, liquid	13390602	31358 [†]	20
Caelo	Wool alcohol ointment SR wit...	14032501	31448	40
Euro OTC	Wool alcohol ointment DAB	1211017-01	31522 [†]	20
Euro OTC	Paraffin, yellow	1401017-01	31541 [†]	20
Caelo	Hydrophilic ointment DAB	14024802	31567 [†]	20
Caelo	Emulsifying hydrophobic base...	14276401	31570 [†]	20
Euro OTC	Paraffin, yellow	1401017-01	31586	60
Euro OTC	Paraffin, liquid	1409007-01	31587 [†]	20
Euro OTC	Paraffin, white	1405014-01	31590 [†]	20
Bombastus	Paraffin, white	290271	31601 [†]	20
Ichthyol	Unguentum Cordes [®]	13A067	31602 [†]	20
Caelo	Eye paraffin	13338711	31610 [†]	20
Caelo	Hydrophilic ointment DAB	13245302	31611 [†]	20
Fagron	Hydrophilic ointment DAB	14G08-T08-013954	31612 [†]	20
Euro OTC	Paraffin, liquid	1506021	31702 [†]	20
Euro OTC	Paraffin, white	1509019	31736 [†]	10
Euro OTC	Paraffin, liquid	1509011	31740 [†]	20
apomix	Wool alcohol ointment SR/DAC	02A070515	31744 [†]	20
Ichthyol	Unguentum Cordes [®]	15A137	31807 [†]	20
Euro OTC	Paraffin, liquid	1511030	31812 [†]	20
Euro OTC	Paraffin, liquid	1511029	31816 [†]	20
Caelo	Hydrophilic ointment DAB	154164	31909 [†]	20
Caelo	Hydrophilic ointment DAB	154165	31910 [†]	20
Caelo	Paraffin, white	153276	31916 [†]	20
Caelo	Paraffin, white	153368	31917 [†]	20
Caelo	Paraffin, white	153368	31918	60
Caelo	Paraffin, white	153497	31919	60
Caelo	Emulsifying hydrophobic base...	153311	31920 [†]	20
Caelo	Wool alcohol ointment SR wit...	153347	31921 [†]	20
Caelo	Paraffin, yellow	153304	31940 [†]	20
Caelo	Paraffin, yellow	153305	31941 [†]	20
Caelo	Paraffin, yellow	153306	31942	60
apomix	Wool alcohol ointment SR/DAC	01A180215	31949 [†]	20
Ichthyol	Unguentum Cordes [®]	15A150	31954 [†]	20
Caelo	Wool alcohol ointment DAB	153423	32011	60
Caelo	Wool alcohol ointment DAB	153424	32012	60
Caelo	Wool alcohol ointment DAB	153427	32013	60
Caelo	Wool alcohol ointment DAB	153428	32014	60
Caelo	Wool alcohol ointment DAB	153429	32015	60
Caelo	Wool alcohol ointment DAB	153558	32016	60
Caelo	Wool alcohol ointment DAB	153562	32017 [†]	20
Caelo	Wool alcohol ointment DAB	153563	32018	60
Caelo	Wool alcohol ointment DAB	153619	32019	60
Caelo	Wool alcohol ointment DAB	153618	32020	60
Caelo	Wool alcohol ointment DAB	153620	32021	48
Caelo	Wool alcohol ointment DAB	161102	32098	40
Caelo	Wool alcohol ointment DAB	161101	32099	40
Caelo	Wool alcohol ointment DAB	161103	32100	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Wool alcohol ointment DAB	161104	32101	40
Caelo	Wool alcohol ointment DAB	161105	32102	40
Caelo	Wool alcohol ointment DAB	161106	32103	40
Caelo	Wool alcohol ointment DAB	161107	32104	40
Euro OTC	Paraffin, liquid	1602026	32189	40
Euro OTC	Paraffin, white	1602009	32191	40
Caelo	Paraffin, white	161380	32223	40
Caelo	Paraffin, white	161380	32224	40
Euro OTC	Paraffin, white	1605019	32282	40
Euro OTC	Paraffin, liquid	1604045	32285	40
Caelo	Paraffin, liquid	161710	32323	40
Caelo	Hydrophilic ointment DAB	161668	32389	40
Caelo	Hydrophilic ointment DAB	161668	32390	40
Caelo	Hydrophilic ointment DAB	161670	32391	40
Euro OTC	Paraffin, liquid	1606032	32462	40
Caelo	Wool alcohol ointment DAB	162488	32511	40
Caelo	Wool alcohol ointment DAB	162488	32512	40
Caelo	Emulsifying hydrophobic base...	162259	32518	40
Caelo	Emulsifying hydrophobic base...	162260	32519	40
Euro OTC	Paraffin, liquid	1608007	32706	40
Euro OTC	Paraffin, white	1608041	32711	40
Euro OTC	Wool alcohol ointment DAB	1608040	32777	40
Euro OTC	Paraffin, liquid	1610017	32830	40
Caelo	Wool alcohol ointment DAB	162951	32864	40
Caelo	Wool alcohol ointment DAB	162952	32865	40
Caelo	Hydrophilic ointment DAB	162824	32866	40
Caelo	Hydrophilic ointment DAB	162825	32867	40
Caelo	Hydrophilic ointment DAB	162826	32868	40
Caelo	Hydrophilic ointment DAB	162827	32869	40
Caelo	Hydrophilic ointment DAB	162828	32870	40
Caelo	Paraffin, white	162411	32901	40
Ichthyol	Unguentum Cordes®	16A061	32912	40
Euro OTC	Paraffin, white	1612007	32982	40
Euro OTC	Paraffin, liquid	1701052	32985	40
Euro OTC	Paraffin, yellow	1701027	32987	40
Caelo	Emulsifying hydrophobic base...	170321	33124	40
Caelo	Emulsifying hydrophobic base...	170322	33132	40
Caelo	Hydrophilic ointment DAB	170403	33134	40
Caelo	Hydrophilic ointment DAB	170406	33135	40
Caelo	Hydrophilic ointment DAB	170402	33159	40
Caelo	Hydrophilic ointment DAB	170405	33160	40
Caelo	Wool alcohol ointment DAB	163318	33184	40
Caelo	Wool alcohol ointment DAB	163317	33185	40
Caelo	Hydrophobic base gel DAC	16151105	33223	40
Euro OTC	Paraffin, liquid	1701033	33236	40
Caelo	Emulsifying hydrophobic base...	171253	33398	40
Caelo	Emulsifying hydrophobic base...	171137	33399	40
Euro OTC	Paraffin, yellow	1705016	33444	42
Euro OTC	Paraffin, white	1705015	33445	40
Caelo	Eye paraffin	15441514	33491	40
apomix	Eye ointment, emulsifying DAC	08A150517	33503	40
apomix	Eye ointment, emulsifying DAC	13C071216	33504	40
Caelo	Emulsifying hydrophobic base...	171979	33545	40
Caelo	Wool alcohol ointment DAB	172062	33547	40
Caelo	Wool alcohol ointment DAB	171416	33548	40
Caelo	Hydrophilic ointment DAB	171855	33554	40
Caelo	Hydrophilic ointment DAB	171563	33555	40
Caelo	Hydrophilic ointment DAB	171564	33556	40
Euro OTC	Paraffin, liquid	1706017	33598	40
Euro OTC	Paraffin, liquid	1707002	33599	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Wool alcohol ointment SR wit...	16251301	33614	40
apomix	Wool alcohol ointment SR/DAC	01B310517	33615	40
apomix	Wool alcohol ointment SR/DAC	01B310517	33616	40
Euro OTC	Paraffin, white	1709016	33659	40
Euro OTC	Paraffin, liquid	1709002	33662	40
Euro OTC	Paraffin, yellow	1709015	33663	40
Euro OTC	Paraffin, liquid	1710025	33848	40
Euro OTC	Paraffin, white	1712008	33934	40
Euro OTC	Paraffin, liquid	1801041	33998	40
Euro OTC	Paraffin, liquid	1802016	34003	40
Caelo	Hydrophobic base gel DAC	173204	34018	40
Caelo	Hydrophilic ointment DAB	170830	34025	40
Caelo	Hydrophilic ointment DAB	172906	34026	40
Caelo	Hydrophilic ointment DAB	173302	34027	40
Euro OTC	Paraffin, white	1804011	34098	40
Euro OTC	Paraffin, yellow	1804012	34099	40
Caelo	Wool alcohol ointment SR wit...	180745	34143	40
Caelo	Emulsifying hydrophobic base...	181406	34155	40
Caelo	Paraffin, yellow	180920	34161	40
Euro OTC	Paraffin, liquid	1806004	34192	40
Caelo	Paraffin, yellow	181915	34255	40
Caelo	Hydrophilic ointment DAB	182075	34261	40
Caelo	Wool alcohol ointment SR wit...	181899	34271	40

- 88265 spectra from a total of 2024 batches from further 425 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 9560 spectra from 464 *Apo-Ident* customers from 2536 batches from the substance/substance group *Paraffin classification*.
- These include spectra of independent samples from 2486 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Apomix/Kehr	Eye ointment, emulsifying DAC	1B020117	1
PKH GmbH Halle	Eye ointment, emulsifying DAC	1C080114	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	6D090516	1
Apomix	Eye ointment, emulsifying DAC	2C040116	1
Apomix/Kehr	Eye ointment, emulsifying DAC	10A041016	1
Apomix/Kehr	Eye ointment, emulsifying DAC	10B041016	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	13A071216	1
Apomix / Phönix	Eye ointment, emulsifying DAC	13B071216	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Eye ointment, emulsifying DAC	14A011013	1
Apomix / Noweda	Eye ointment, emulsifying DAC	15A030815	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	16A231014	1
AMH	Eye ointment, emulsifying DAC	16B011015	1
PHOENIX 03.03.2014	Eye ointment, emulsifying DAC	16B111113	1
Apomix/Noweda	Eye ointment, emulsifying DAC	17A231014	1
PhR Neidhardt	Eye ointment, emulsifying DAC	20A011215	1
Noweda/apomix	Eye ointment, emulsifying DAC	1B070113	1
apomix/Noweda	Eye ointment, emulsifying DAC	20D011215	1
Caelo	Eye ointment, emulsifying DAC	3A180213	1
Apomix / Phönix	Eye ointment, emulsifying DAC	7A210513	1
PKH Halle / Noweda	Eye ointment, emulsifying DAC	7D210513	1
apomix/Phönix	Eye ointment, emulsifying DAC	9B25061	1
Apomix	Eye ointment, emulsifying DAC	11B030912	1
Caelo	Eye ointment, emulsifying DAC	13084905	1
Alliance Healthcare/Apo...	Eye ointment, emulsifying DAC	17A111113	1
Caelo	Eye paraffin	5011507	1
Caelo	Eye paraffin	6021505	1
apomix	Eye ointment, emulsifying DAC	17B111113	1
Caelo	Eye paraffin	13338710	1
Caelo	Eye paraffin	13338707	1
Gehe	Eye paraffin	13338704	1
Caelo	Eye paraffin	13338711	3
Caelo	Eye paraffin	13338713	10
Caelo	Eye paraffin	14124703	6
Caelo	Eye paraffin	14124704	6
Sanacorp	Eye paraffin	14124704	1
Phoenix Goettin	Eye paraffin	14124704	1
Caelo	Eye paraffin	14124705	7
Sanacorp	Eye paraffin	14124705	1
Caelo	Eye paraffin	14124706	6
Gehe	Eye paraffin	14124706	1
Caelo	Eye paraffin	14124707	16
Caelo	Eye paraffin	14124708	5
Caelo	Eye paraffin	1412704	1
Caelo	Eye paraffin	1412703	1
Caelo	Eye paraffin	14301804	12
Caelo	Eye paraffin	14301805	10
Caelo/ Noweda	Eye paraffin	14301804	1
Caesar & Loretz GmbH / ...	Eye paraffin	14301804	1
Caelo	Eye paraffin	15050002	27
Caelo	Eye paraffin	14551510	1
Caelo	Eye paraffin	14301806	1
Caesar & Loretz/Noweda	Eye paraffin	150500004	1
Caele	Eye paraffin	15050002	1
Caelo	Eye paraffin	15050003	12
Caelo	Eye paraffin	15050004	6
Caelo	Eye paraffin	15050005	6
Caelo	Eye paraffin	15050007	3
Caelo	Eye paraffin	1505002	2
Caelo	Eye paraffin	1505003	1
Caelo	Eye paraffin	1527880	1
Caelo	Eye paraffin	152778809	1
Caelo	Eye paraffin	15278801	15
Caelo	Eye paraffin	15278802	2
Caelo	Eye paraffin	15278803	2
Caelo	Eye paraffin	15278804	13
Caelo	Eye paraffin	15278803/170314VA	1
Caelo	Eye paraffin	15278805	12
Caelo	Eye paraffin	15278808	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Eye paraffin	15278809	5
Ebert & Jacobi 31.05...	Eye paraffin	15278809	1
Caelo	Eye paraffin	1544111	1
Caelo	Eye paraffin	15441414	1
Caelo	Eye paraffin	15441501	8
Caelo	Eye paraffin	15441502	19
Caelo	Eye paraffin	15441504	18
Caecolo/Noweda	Eye paraffin	15441504	1
Caelo	Eye paraffin	15441505	18
Ebert u. Jacobi 20.07...	Eye paraffin	15441505	1
Caelo	Eye paraffin		3
Caelo	Eye paraffin	15441506	4
Caelo	Eye paraffin	15441509	16
Caelo	Eye paraffin	15441510	6
Noweda	Eye paraffin	15441510	2
Caelo	Eye paraffin	15441511	15
Gehe	Eye paraffin	15441511	1
Caelo	Eye paraffin	15441514	19
phönix	Eye paraffin	15441514	1
Fagron	Eye paraffin	15441514	2
Caecolo	Eye paraffin	15441514	1
Euro OTC	Eye paraffin	15441514	1
Caelo	Eye paraffin	15441515	1
Caelo	Eye paraffin	1544514	2
Caelo	Eye paraffin	160427VA0	1
Caelo	Eye paraffin	16096913	2
Caelo	Eye paraffin	16121502	1
Finteler Apotheke	Eye paraffin	16271807	2
Caelo	Eye paraffin	16271807	7
Noweda	Eye paraffin	16271807	1
Caelo	Eye paraffin	16271808	2
Caelo	Eye paraffin	16271810	7
Caelo	Eye paraffin	17052001	7
Caelo	Eye paraffin	17052003	1
Caelo	Eye paraffin	17052008	2
Caelo	Eye paraffin	171123VPA/16271808	1
Caelo	Eye paraffin	19/AKD13	1
Caelo	Eye paraffin	180206AV/17052002	1
Caelo	Eye paraffin	19/AKD2	1
Sanacorp	Eye paraffin	19/AKD4	1
Caelo	Eye paraffin	19/AKD7	1
Caelo	Eye paraffin	24031503	1
Caelo	Eye paraffin	29101401	1
Caelo	Eye paraffin	19051504	1
Caelo	Eye paraffin	20101401	1
Caelo	Eye paraffin	31071502	1
Fagron	Eye paraffin	760619-021159	1
Fagron	Eye paraffin	790619-013771	1
Fagron	Eye paraffin	790619-021159	3
Sanacorp	Eye paraffin		1
Bombastus	Hydrophilic ointment DAB	284748	5
Caelo	Paraffin, white		5
gehe	Hydrophobic base gel DAC	10L14-T54-006109	3
Caelo	Hydrophobic base gel DAC	13341703	3
Caelo	Hydrophobic base gel DAC	14044601	12
Caelo	Hydrophobic base gel DAC	14087501	1
Caelo	Hydrophobic base gel DAC	14102502	5
Caelo	Hydrophobic base gel DAC	14158002	4
Caelo	Hydrophobic base gel DAC	14158003	9
Caelo	Hydrophobic base gel DAC	14158006	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophobic base gel DAC	14161501	1
Caelo	Hydrophobic base gel DAC	14276401	3
Caelo	Hydrophobic base gel DAC	15055705	5
Caelo	Hydrophobic base gel DAC	15064101	8
Caelo	Hydrophobic base gel DAC	15149603	3
Caelo	Hydrophobic base gel DAC	15064201	1
Caelo	Hydrophobic base gel DAC	15205002	1
Caelo	Hydrophobic base gel DAC	15149604	8
Caelo	Hydrophobic base gel DAC	15149605	3
Caelo	Hydrophobic base gel DAC	15240801	7
Caelo	Hydrophobic base gel DAC	15240805	7
Caelo	Hydrophobic base gel DAC	15269201	1
Caelo	Hydrophobic base gel DAC	15331102	1
Caelo	Hydrophobic base gel DAC	15436103	2
Caelo	Hydrophobic base gel DAC	15436106	1
Caelo	Hydrophobic base gel DAC	16054201	5
Caelo	Hydrophobic base gel DAC	16070501	1
Caelo	Hydrophobic base gel DAC	16070502	5
Caelo	Hydrophobic base gel DAC	16070503	9
Caesar & Loretz GmbH Hi...	Hydrophobic base gel DAC	16070503	2
Caelo	Hydrophobic base gel DAC	16070505	6
Caelo	Hydrophobic base gel DAC	16151107	1
Caelo	Hydrophobic base gel DAC	16172401	1
Caelo	Hydrophobic base gel DAC	16240104	4
Fagron	Hydrophobic base gel DAC	16H24-T02-034643	3
Fagron	Hydrophobic base gel DAC	17032201	1
Bombastus	Hydrophobic base gel DAC	18001099	1
Caesar&Loretz GmbH	Hydrophobic base gel DAC	17043605	1
Bombastus	Hydrophobic base gel DAC	28111403	1
Bombastus	Hydrophobic base gel DAC	282094	5
Caelo	Hydrophobic base gel DAC	282094	2
GEHE 02.12.2014	Hydrophobic base gel DAC	286013	2
Bombastus	Hydrophobic base gel DAC	289213	1
Bombastus	Hydrophobic base gel DAC	294417	1
Bombastus	Hydrophobic base gel DAC	289214	5
Bombastus	Hydrophobic base gel DAC	294418	1
Bombastus	Hydrophobic base gel DAC	295729	4
Bombastus	Hydrophobic base gel DAC	300896	2
Bombastus	Hydrophobic base gel DAC	304036	4
Fagron	Hydrophobic base gel DAC	10L14-T54	3
Caelo	Hydrophobic base gel DAC	9082013C	1
Fagron	Hydrophobic base gel DAC	10L14-T54-000713	1
Fagron	Hydrophobic base gel DAC	11K23-T01	1
Caelo	Hydrophobic base gel DAC	12229001	2
Caelo	Hydrophobic base gel DAC	12229003	4
Caelo	Hydrophobic base gel DAC	12229005	1
Caelo	Hydrophobic base gel DAC	12362503	5
11.07.13	Hydrophobic base gel DAC	12324704	1
Caelo	Hydrophobic base gel DAC	12362504	2
Caelo	Hydrophobic base gel DAC	12362505	3
Euro OTC	Hydrophobic base gel DAC	12362505	1
Caelo	Hydrophobic base gel DAC	13085004	3
Caelo	Hydrophobic base gel DAC	123662505	1
Caesar & Loretz GmbH, 4...	Hydrophobic base gel DAC	13165201	1
Caelo	Hydrophobic base gel DAC	13165201	3
Europ Ots Pharma	Hydrophobic base gel DAC	13299604	1
Caelo	Hydrophobic base gel DAC	13299604	10
Caelo	Hydrophobic base gel DAC	13295801	1
Caelo	Hydrophobic base gel DAC	14102504	2
15.07.13	Hydrophobic base gel DAC	277443	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrophobic base gel DAC	277443	2
Sanacorp	Hydrophobic base gel DAC	AKD18	1
Bombastus	Hydrophobic base gel DAC	279074	4
Caelo	Emulsifying hydrophobic base...	30071201	1
Ichthyol	Unguentum Cordes®	4091516	1
Ichthyol	Unguentum Cordes®	5052014B	1
Ichthyol	Unguentum Cordes®	7011501	1
Sanacorp	Unguentum Cordes®	103/AKD1	1
Ichthyol /Noweda	Unguentum Cordes®	12A118	1
Caelo	Unguentum Cordes®	13A042	2
Ichthyol	Unguentum Cordes®	13012015B	1
Ichthyol /Noweda	Unguentum Cordes®	13A067	1
Ichthyol	Unguentum Cordes®	13041504	1
Ichthyol / Gehe	Unguentum Cordes®	13A068	2
ichthyol-gesellschaft	Unguentum Cordes®	13A068	2
Ichthyol/Noweda	Unguentum Cordes®	13A089	2
Ichthyol/Phönix	Unguentum Cordes®	13A089	1
Alliance Health 11.03.14	Unguentum Cordes®	13A090	1
Caesar & Loretz GmbH Hi...	Unguentum Cordes®	13A090	1
Caelo	Unguentum Cordes®	13A090	2
Ichthyol/Noweda	Unguentum Cordes®	13A090	2
Ichthyol / Phoenix	Unguentum Cordes®	13A090	1
ichthyol/Sanacorp	Unguentum Cordes®	13A090	1
Noweda	Unguentum Cordes®	13A091	1
IchthyolSanacorp	Unguentum Cordes®	13A110	1
Ichthyol	Unguentum Cordes®	13A138	3
Caelo	Unguentum Cordes®	13A136	1
ichthyol/Gehe	Unguentum Cordes®	13A091	1
Ichthyol	Unguentum Cordes®	13A140	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	13A138	1
Ichthyol/Phönix	Unguentum Cordes®	13A140	1
Phoenix/Ichthyol	Unguentum Cordes®	13A140	1
Ichthyol	Unguentum Cordes®	13A141	1
Caelo	Unguentum Cordes®	13A141	1
Ichthyol/Sanacorp	Unguentum Cordes®	13A142	1
Ichthyol / Gehe	Unguentum Cordes®	13A142	3
Ichthyol-Gesellschaft	Unguentum Cordes®	13A142	1
Ichthyol / Gehe	Unguentum Cordes®	13A153	1
Ichthyol	Unguentum Cordes®	13A142	1
Ichthyol/GEHE	Unguentum Cordes®	13A142	1
Caelo	Unguentum Cordes®	13A153	2
Ichthyol	Unguentum Cordes®	13A153	1
Caelo	Unguentum Cordes®	13A157	1
Ichthyol / Gehe	Unguentum Cordes®	13A156	1
Ichthyol Gesellschaft	Unguentum Cordes®	13A156	1
Ichthyol/Noweda	Unguentum Cordes®	13A156	1
Gehe	Unguentum Cordes®	13A156	1
Cordes/Gehe	Unguentum Cordes®	13A157	1
Ichthyol/Noweda	Unguentum Cordes®	13A158	1
Berkel AHK	Unguentum Cordes®	13A158	1
Caelo	Unguentum Cordes®	13A159	1
Ichthyol/Ebert	Unguentum Cordes®	13A159	1
Sanacorp	Unguentum Cordes®	13A159	1
Gehe	Unguentum Cordes®	13A159	1
Ichthyol	Unguentum Cordes®	13A159	2
Ichthyol/Otto	Unguentum Cordes®	13A159	1
Ichthyol/Phoenix	Unguentum Cordes®	13A166	1
Ichthyol/Phönix	Unguentum Cordes®	13A166	3
Ichthyol / Gehe	Unguentum Cordes®	13A166	1
Ichthyol	Unguentum Cordes®	13A166	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanova Pharma	Unguentum Cordes [®]	13A166	1
Ichthyol/Fiebig	Unguentum Cordes [®]	13A166	1
Ichthyol über Phoenix	Unguentum Cordes [®]	13A166	1
Caelo	Unguentum Cordes [®]	13A171	2
Ichthyol Gesellschaft	Unguentum Cordes [®]	13A170	1
Ichthyol	Unguentum Cordes [®]	13A171	1
Ichthyol	Unguentum Cordes [®]	13A172	2
Anzag	Unguentum Cordes [®]	13A172	1
Ichthyol/Noweda	Unguentum Cordes [®]	13A172	1
Cordes/Gehe	Unguentum Cordes [®]	13A172	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	13A172	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	13A172	1
Caelo	Unguentum Cordes [®]	13A172	1
Ichthyol	Unguentum Cordes [®]	13A178	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	13A178	1
Ichthyol/Gehe	Unguentum Cordes [®]	13A178	3
Ichthyol-Gesellschaft	Unguentum Cordes [®]	13A180	3
Anzag/Ichthyol	Unguentum Cordes [®]	13A180	1
Ichthyol	Unguentum Cordes [®]	13A180	1
Ichthyol / Gehe	Unguentum Cordes [®]	13A181	1
Ichthyol -Gesellschaft	Unguentum Cordes [®]	13A181	1
Ichthyol/Gehe	Unguentum Cordes [®]	13A181	2
Caelo	Unguentum Cordes [®]	13A181	2
Ichthyol Gesellschaft/ ...	Unguentum Cordes [®]	13A181	1
Ichthyol /Noweda	Unguentum Cordes [®]	13A181	1
Ichthyol/Noweda	Unguentum Cordes [®]	13A181	1
Ichthyol/Gehe	Unguentum Cordes [®]	13A216	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	13A214	1
ANZAG/Ichthyol Gesellsc...	Unguentum Cordes [®]	13A72	1
Gehe/Ichthyol	Unguentum Cordes [®]	13A090	1
Ichthyol	Unguentum Cordes [®]	14031409	1
Ichthyol/Anzag	Unguentum Cordes [®]	1415A-02532	1
Ichthyol-Gesellschaft/ ...	Unguentum Cordes [®]	14A005	1
Ichthyol/Gehe	Unguentum Cordes [®]	14A006	1
Ichthyol Cordes/Allianc...	Unguentum Cordes [®]	14A009	1
Caelo	Unguentum Cordes [®]	14A009	1
Ichthyol/Noweda	Unguentum Cordes [®]	14A009	1
Ichthyol/Noweda	Unguentum Cordes [®]	14A047	1
Ichthyol Gesellschaft/ ...	Unguentum Cordes [®]	14A047	1
Ichthyol/Gehe	Unguentum Cordes [®]	14A055	2
Caelo	Unguentum Cordes [®]	14A055	2
Gehe	Unguentum Cordes [®]	14A056	1
Ichthyol/ Gehe	Unguentum Cordes [®]	14A056	1
Ichthyol/Gehe	Unguentum Cordes [®]	14A056	2
Cordes	Unguentum Cordes [®]	14A055	1
Berkel AHK	Unguentum Cordes [®]	14A056	1
Stumpf/Ichthyol	Unguentum Cordes [®]	14A057	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	14A057	1
Ichthyol	Unguentum Cordes [®]	14A057	3
Ichthyol/Phönix	Unguentum Cordes [®]	14A057	1
sanacorp	Unguentum Cordes [®]	14A060	1
Bombastus	Unguentum Cordes [®]	14A060	1
Ichthyol	Unguentum Cordes [®]	14A060	1
Ichthyol/Ebert	Unguentum Cordes [®]	14A060	1
sanacorp/ichthyol	Unguentum Cordes [®]	14A061	2
Ichthyol	Unguentum Cordes [®]	14A061	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	14A061	2
Gehe/Ichthyol	Unguentum Cordes [®]	14A061	1
Ichthyol/Gehe	Unguentum Cordes [®]	14A062	3
Cordes	Unguentum Cordes [®]	14A062	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol Gesellschaft C...	Unguentum Cordes®	14A062	1
Caelo	Unguentum Cordes®	14A063	2
Euro OTC	Unguentum Cordes®	14A078	2
Ichthyol-Gesellschaft	Unguentum Cordes®	14A063	1
Ichthyol	Unguentum Cordes®	14A063	1
Stumpf/Ichtyol	Unguentum Cordes®	14A063	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	14A063	2
Ichthyol-Gesellschaft C...	Unguentum Cordes®	14A063	1
Ichthyol über Phoenix	Unguentum Cordes®	14A062	1
Ichthyol/Phönix	Unguentum Cordes®	14A062	1
anzag	Unguentum Cordes®	14A079	1
Caelo	Unguentum Cordes®	14A078	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	14A079	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	14A079	2
Euro OTC	Unguentum Cordes®	14A079	2
Ichthyol / E	Unguentum Cordes®	14A079	1
Alliance health	Unguentum Cordes®	14A079	1
Ichthyol/Phoenix	Unguentum Cordes®	14A079	1
Ichthyol cordes / Noweda	Unguentum Cordes®	14A079	1
Ichtyol-Gesellschaft	Unguentum Cordes®	14A079	1
eurho	Unguentum Cordes®	14a079	1
Ichthyol-Gesellschaft	Unguentum Cordes®	14A079	1
Caelo	Unguentum Cordes®	14A079	1
Ichthyol-Gesellschaft/P...	Unguentum Cordes®	14A095	2
Ichtyol-Gesellschaft/Ph...	Unguentum Cordes®	14A095	1
Ichthyol/Geilenkirchen	Unguentum Cordes®	14A095	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	14A095	1
Ichthyol-Gesellschaft	Unguentum Cordes®	14A095	1
Ichthyol/Phönix	Unguentum Cordes®	14A095	1
Ichthyol	Unguentum Cordes®	14A095	1
Ichtyol/Sanacorp	Unguentum Cordes®	14A097	1
Ichthyol	Unguentum Cordes®	14A098	1
Ichthyol/ Noweda	Unguentum Cordes®	14A099	1
Ichthyol/Sanacorp	Unguentum Cordes®	14A098	1
Ichtyol - Gehe	Unguentum Cordes®	14A1115	1
Ichthyol/Gehe	Unguentum Cordes®	14A112	1
Ichthyol	Unguentum Cordes®	14A113	1
Caelo	Unguentum Cordes®	14A113	2
Noweda/Ichthyol	Unguentum Cordes®	14A113	1
Ichthyol/Otto	Unguentum Cordes®	14A113	1
Cealo	Unguentum Cordes®	14A113	1
Ichthyol/GEHE	Unguentum Cordes®	14a113	1
Ichthyol-Gesellschaft	Unguentum Cordes®	14A114	2
Ichthyol	Unguentum Cordes®	14A114	1
Ichthyol über Phoenix	Unguentum Cordes®	14A114	1
ICHTHYOL	Unguentum Cordes®	14A114	1
Ichthyol-Gesellschaft/P...	Unguentum Cordes®	14a114	1
Ichthyol-Gesellschaft	Unguentum Cordes®	14A115	2
Caelo	Unguentum Cordes®	14A116	1
Ichthyol/Phönix	Unguentum Cordes®	14A116	1
Cordes	Unguentum Cordes®	14A117	2
phönix	Unguentum Cordes®	14A116	1
AHD	Unguentum Cordes®	14A117	1
sanacorp/ichthyol	Unguentum Cordes®	14A117	1
Ichthyol/Gehe	Unguentum Cordes®	14A117	1
Ichthyol/Gehe	Unguentum Cordes®	14A130	1
Ichthyol/Sanacorp	Unguentum Cordes®	14A171	1
Ichthyol / Sanacorp	Unguentum Cordes®	14A168	1
Ichthyol	Unguentum Cordes®	14A132	1
ICHTHYOL	Unguentum Cordes®	14A171	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol / Noweda	Unguentum Cordes [®]	14A172	1
Ichthyol	Unguentum Cordes [®]	14A171	2
Caelo	Unguentum Cordes [®]	14A172	1
Ichthyol / Sanacorp	Unguentum Cordes [®]	14A172	1
Ichthyol	Unguentum Cordes [®]	14A172	2
Sanacorp	Unguentum Cordes [®]	14A172	1
Ichthyol- Gesellschaft ...	Unguentum Cordes [®]	14A172	1
Ichthyol/Phönix	Unguentum Cordes [®]	14A172	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	14A172	2
Caelo	Unguentum Cordes [®]	14A173	1
Ichthyol/ Noweda	Unguentum Cordes [®]	14A172	1
ichthyol	Unguentum Cordes [®]	14A172	1
Ichthyol / Noweda	Unguentum Cordes [®]	14A173	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	14A173	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	14A179	1
Ichthyol	Unguentum Cordes [®]	14A173	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	14A179	1
Ichthyol/Fiebig	Unguentum Cordes [®]	14A179	1
Ichthyol/Phönix	Unguentum Cordes [®]	14A179	1
Caelo	Unguentum Cordes [®]	14A213	2
Förster/Phönix	Unguentum Cordes [®]	158142	1
Caelo	Unguentum Cordes [®]	14A312	1
Fagron	Unguentum Cordes [®]	15A010	1
Ichthyol/Otto	Unguentum Cordes [®]	15A010	1
ICHTHYOL/GEHE	Unguentum Cordes [®]	15A011	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A011	1
Ichthyol Pharma	Unguentum Cordes [®]	15A011	1
Ichthyol-Gesellschaft /...	Unguentum Cordes [®]	15A020	2
Ichthyol/Anzag	Unguentum Cordes [®]	15A020	1
ICHTHYOL Gesellschaft /...	Unguentum Cordes [®]	15A020	1
Fagron	Unguentum Cordes [®]	15A020	1
Bombastus	Unguentum Cordes [®]	15a020	1
Ichthyol/phönix	Unguentum Cordes [®]	15A020	1
Ichtyol-Gesellschaft	Unguentum Cordes [®]	15A020	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A020	2
Ichthyol/ Gehe	Unguentum Cordes [®]	15A020	2
Ichthyol / G	Unguentum Cordes [®]	15A020	1
Caelo	Unguentum Cordes [®]	15A020	1
Caelo	Unguentum Cordes [®]	15a020	1
Anzag	Unguentum Cordes [®]	15A021	1
Ichthyol / Noweda	Unguentum Cordes [®]	15A021	1
Gehe Kassel	Unguentum Cordes [®]	15A021	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A022	1
Ichthyol/ Noweda	Unguentum Cordes [®]	15a021	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A022	1
sanacorp/ichthyol	Unguentum Cordes [®]	15A022	1
Anzag	Unguentum Cordes [®]	15A023	1
Ichthyol	Unguentum Cordes [®]	15A022	1
Ichthyol-Gesellschaft /...	Unguentum Cordes [®]	15A023	1
Ichthyol, Noweda	Unguentum Cordes [®]	15A023	1
Caelo	Unguentum Cordes [®]	15A027	1
Gehe - Ichthyol	Unguentum Cordes [®]	15A027	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A027	1
Ichthyol	Unguentum Cordes [®]	15A023	1
Caelo	Unguentum Cordes [®]	15A023	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A027	4
Caelo	Unguentum Cordes [®]	15A040	1
Gehe,11.05.15/Cordes	Unguentum Cordes [®]	15A040	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A040	2
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A037	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/ Noweda	Unguentum Cordes [®]	15A040	1
Ichthyol/Fiebig	Unguentum Cordes [®]	15A040	1
Ichthyol	Unguentum Cordes [®]	15A040	3
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A041	1
Ichthyol GmbH	Unguentum Cordes [®]	15A041	1
Ichthyol über Phoenix	Unguentum Cordes [®]	15A042	1
Ichthyol Pharma	Unguentum Cordes [®]	15A042	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A042	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A042	1
sanacorp/ichthyol	Unguentum Cordes [®]	15A042	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A042	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A043	2
Sanacorp	Unguentum Cordes [®]	15A043	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	15A043	1
sanacorp/ichthyol	Unguentum Cordes [®]	15A044	1
SANOVA	Unguentum Cordes [®]	15A044	1
Caelo	Unguentum Cordes [®]	15A060	1
Ichthyol-Gesellschaft /...	Unguentum Cordes [®]	15A061	1
Ichthyol /Noweda	Unguentum Cordes [®]	15A060	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	15A061	1
Noweda	Unguentum Cordes [®]	15A060	1
ichthyol	Unguentum Cordes [®]	15A061	1
anzag	Unguentum Cordes [®]	15A062	1
Caelo	Unguentum Cordes [®]	15A062	2
Noweda	Unguentum Cordes [®]	15A062	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A062	1
Ichthyol / Noweda	Unguentum Cordes [®]	15A073	3
Ichthyol/Phönix	Unguentum Cordes [®]	15A073	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	15A073	1
Caelo	Unguentum Cordes [®]	15A073	2
Ichthyol Gesellschaft/ ...	Unguentum Cordes [®]	15A087	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A074	1
Ichthtyol	Unguentum Cordes [®]	15A089	1
Ichthyol	Unguentum Cordes [®]	15A079	1
IchthyolSanacorp	Unguentum Cordes [®]	15A076	1
ICHTHYOL	Unguentum Cordes [®]	15a089	1
Caelo	Unguentum Cordes [®]	15A089	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A089	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A090	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A090	1
Ichthyol	Unguentum Cordes [®]	15A091	2
Ichthyol/Phoenix	Unguentum Cordes [®]	15A090	1
Caelo	Unguentum Cordes [®]	15A091	1
Caelo	Unguentum Cordes [®]	15A092	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A091	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A092	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A092	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A106	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A108	3
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A106	1
Ichthyol/Otto	Unguentum Cordes [®]	15A107	1
sanacorp/ichthyol	Unguentum Cordes [®]	15A108	1
Ichthyol	Unguentum Cordes [®]	15A108	1
Noweda	Unguentum Cordes [®]	15A123	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	15A12	1
Ichthyol /Noweda	Unguentum Cordes [®]	15A123	3
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	15A123	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A123	1
Fiebig/Ichthyol	Unguentum Cordes [®]	15A124	1
Ichthyol-Gesellschaft /...	Unguentum Cordes [®]	15A124	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe/Ichthyol	Unguentum Cordes [®]	15A124	4
Ichthyol/Noweda	Unguentum Cordes [®]	15A124	1
ichthyol	Unguentum Cordes [®]	15A124	2
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	15A124	1
Ichthyol/Phoenix WE:02...	Unguentum Cordes [®]	15A137	1
ichthyol	Unguentum Cordes [®]	15A125	1
Caelo	Unguentum Cordes [®]	15A137	2
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A137	1
Sanacorp	Unguentum Cordes [®]	15A137	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A137	1
Ichthyol Pharma	Unguentum Cordes [®]	15A139	1
Ichthyol/ Anzag	Unguentum Cordes [®]	15A139	1
ichthyol	Unguentum Cordes [®]	15A139	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A139	1
Caelo	Unguentum Cordes [®]	15A139	1
VDL	Unguentum Cordes [®]	15A140	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A140	2
Ichthyol/Phönix	Unguentum Cordes [®]	15A140	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A140	1
Ichthyol Pharma	Unguentum Cordes [®]	15A141	1
Ichthyol/ Anzag	Unguentum Cordes [®]	15A141	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A142	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A142	1
Ichthyol	Unguentum Cordes [®]	15A142	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A142	1
Ichthyol	Unguentum Cordes [®]	15A148	4
Ichthyol	Unguentum Cordes [®]	15A150	1
Ichthyol /Noweda	Unguentum Cordes [®]	15A148	1
Ichthyol/ Phönix	Unguentum Cordes [®]	15A150	4
Ichthyol/Gehe	Unguentum Cordes [®]	15A148	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A148	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A150	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A150	1
Caelo	Unguentum Cordes [®]	15A155	1
Ichthyol/ Anzag	Unguentum Cordes [®]	15A155	1
Euro OTC	Unguentum Cordes [®]	15A155	1
Ichthyol Pharma	Unguentum Cordes [®]	15A155	1
sanacorp/ichthyol	Unguentum Cordes [®]	15A156	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A156	2
Ichthyol	Unguentum Cordes [®]	15A156	1
Ichthyol /Noweda	Unguentum Cordes [®]	15A159	1
Ichthyol/Sanova Pharma	Unguentum Cordes [®]	15A156	1
Ichthyol	Unguentum Cordes [®]	15A159	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A164	1
Noweda	Unguentum Cordes [®]	15A191	1
Ichthyol	Unguentum Cordes [®]	15A191	3
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A191	1
Ichthyol /Noweda	Unguentum Cordes [®]	15A191	1
Ichthyol / K	Unguentum Cordes [®]	15A191	1
Gehe/Ichthyol	Unguentum Cordes [®]	15A191	4
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A191	1
Ichthyol/GEHE	Unguentum Cordes [®]	15A089	1
Ichthyol/Anzag	Unguentum Cordes [®]	15A195	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A191	2
Caelo	Unguentum Cordes [®]	15A191	2
Phoenix	Unguentum Cordes [®]	15ao93	1
Ichthyol/GEHE	Unguentum Cordes [®]	15A061	1
Caelo	Unguentum Cordes [®]	16/082/C	1
Ichthyol/Phönix	Unguentum Cordes [®]	16A119	1
Ichthyol Gesellschaft/N...	Unguentum Cordes [®]	16a011	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/Phönix	Unguentum Cordes®	16A142	1
Ichthyol-Gesellschaft	Unguentum Cordes®	16A022	1
SANOVA	Unguentum Cordes®	16A045	1
Ichthyol/Sanacorp	Unguentum Cordes®	16A026	1
Ichthyol-Gesellschaft	Unguentum Cordes®	16A044	1
Ichthyol-Gesellschaft	Unguentum Cordes®	16A046	1
Ichthyol	Unguentum Cordes®	16A022	1
Ichthyol Gesellschaft	Unguentum Cordes®	16a045	1
Ichthyol	Unguentum Cordes®	16A045	1
Ichthyol/ Anzag	Unguentum Cordes®	16A046	2
Ichthyol	Unguentum Cordes®	16A047	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	16A047	1
Ichthyol	Unguentum Cordes®	16A048	2
Sanacorp	Unguentum Cordes®	16A048	2
Ichthyol/Anzag	Unguentum Cordes®	16A048	1
Phönix	Unguentum Cordes®	16A048	1
Caelo	Unguentum Cordes®	16A051	1
Ichthyol über Phoenix	Unguentum Cordes®	16A058	1
Ichthyol	Unguentum Cordes®	16A058	3
Euro OTC	Unguentum Cordes®	16A074	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	16A083	1
Caelo	Unguentum Cordes®	16A077	1
Bombastus	Unguentum Cordes®	16A083	1
AHD	Unguentum Cordes®	16A085	1
Ichthyol /Noweda	Unguentum Cordes®	16A085	1
Ichthyol	Unguentum Cordes®	16A085	2
Ichthyol Gesellschaft C...	Unguentum Cordes®	16A099	1
Caelo	Unguentum Cordes®	16A099	2
Ichthyol-Gesellschaft	Unguentum Cordes®	16A085	1
Ichthyol / Phoenix	Unguentum Cordes®	16A099	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	16A099	1
GEHE	Unguentum Cordes®	16A099	1
	Unguentum Cordes®	16A099	2
Bombastus	Unguentum Cordes®	16A099	1
Ichthyol-Gesellschaft /...	Unguentum Cordes®	16A099	1
Ichthyol-Gesellschaft	Unguentum Cordes®	16A101	2
ichthyol-G.	Unguentum Cordes®	16A099	1
Cordes/Noweda	Unguentum Cordes®	16A101	1
Ichthyol Gesellschaft	Unguentum Cordes®	16A101	2
Ichthyol Gesellschaft/N...	Unguentum Cordes®	16A101	1
Ichthyol-Gesellschaft/K...	Unguentum Cordes®	16A101	1
ichthyol	Unguentum Cordes®	16A101	3
Ichthyol-Gesellschaft	Unguentum Cordes®	16A102	2
Ichthyol Gesellschaft	Unguentum Cordes®	16A102	1
Caelo	Unguentum Cordes®	16A102	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	16A102	1
Ichthyol Gesellschaft	Unguentum Cordes®	16A119	1
Ichthyol Gesellschaft/ ...	Unguentum Cordes®	16A102	1
Ichthyol/ Anzag	Unguentum Cordes®	16A119	1
Ichthyol über Phoenix	Unguentum Cordes®	16A119	1
Ichthyol	Unguentum Cordes®	16A119	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	16A121	1
Ichthyol-Gesellschaft	Unguentum Cordes®	16A121	1
Phoenix	Unguentum Cordes®	16A121	1
Ichthyol/Phoenix	Unguentum Cordes®	16A121	4
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes®	16A121	1
Cordes	Unguentum Cordes®	16A122	1
Caelo	Unguentum Cordes®	16A122	2
Ichthyol	Unguentum Cordes®	16A122	1
Ichthyol / Phoenix	Unguentum Cordes®	16A139	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol Gesellschaft	Unguentum Cordes [®]	16A124	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	16A123	1
Ichthyol/ Noweda	Unguentum Cordes [®]	16A137	1
Caelo	Unguentum Cordes [®]	16A139	2
Ichthyol Gesellschaft	Unguentum Cordes [®]	16A139	1
Caelo	Unguentum Cordes [®]	16A140	2
Ichthyol Gesellschaft	Unguentum Cordes [®]	16A140	1
Ichthyol/Phoenix	Unguentum Cordes [®]	16A140	3
Ichthyol/Noweda 31.01.2017	Unguentum Cordes [®]	16A140	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	16A140	1
PhoenixWE:7.02.2017 E...	Unguentum Cordes [®]	16A141	1
Ichthyol/ Anzag	Unguentum Cordes [®]	16A141	1
Ichthyol/Noweda	Unguentum Cordes [®]	16A142	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A142	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	16A142	1
Cordes/Phönix	Unguentum Cordes [®]	16A142	2
Caelo	Unguentum Cordes [®]	16A142	1
Ichthyol	Unguentum Cordes [®]	16A142	1
Ichthyol Gesellschaft/ ...	Unguentum Cordes [®]	16A144	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A143	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A144	1
Caelo	Unguentum Cordes [®]	16A144	2
Cordes/Gehe	Unguentum Cordes [®]	16A155	1
Ichthyol/GEHE	Unguentum Cordes [®]	16A143	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A153	1
Ichthyol/ Phoenix	Unguentum Cordes [®]	16A155	1
Ichthyol/Noweda	Unguentum Cordes [®]	16A15501062983	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	16A173	1
Ichthyol	Unguentum Cordes [®]	16A173	3
Ichthyol	Unguentum Cordes [®]	16A176	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	16A183	1
Ichthyol / G	Unguentum Cordes [®]	16A184	1
Ichthyol-Gesellschaft/K...	Unguentum Cordes [®]	16A184	1
Allianc. Health. / Icht...	Unguentum Cordes [®]	16A184	1
Caelo	Unguentum Cordes [®]	16A184	2
Ichthyol Gesellschaft/P...	Unguentum Cordes [®]	16A184	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A184	2
Caelo	Unguentum Cordes [®]	16A210	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A210	1
Fagron	Unguentum Cordes [®]	1744	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A002	1
Ichthyol/Phönix	Unguentum Cordes [®]	16A048	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	17A002	1
Ichthyol / AHD	Unguentum Cordes [®]	17A002	1
ichthyol	Unguentum Cordes [®]	17A003	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	17A003	1
Ichthyol	Unguentum Cordes [®]	17A053	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A069	1
Fagron	Unguentum Cordes [®]	17a054	1
Caelo	Unguentum Cordes [®]	17A163	1
Caelo	Unguentum Cordes [®]	17A070	1
Ichthyol/Gehe	Unguentum Cordes [®]	17A164	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A164	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A101	1
Caelo	Unguentum Cordes [®]	17A217	1
ichthyol	Unguentum Cordes [®]	17A163	1
Ichthyol	Unguentum Cordes [®]	30092014C	1
Ichthyol	Unguentum Cordes [®]	414022-2-P	1
Ichthyol/Gehe	Unguentum Cordes [®]	4913Q-02532	1
Beiersdorf AG	Unguentum Cordes [®]	62932210	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix /Ichthyol	Unguentum Cordes [®]	Ch13A089	1
Ichtyol	Unguentum Cordes [®]	11A038	1
gebepharm	Unguentum Cordes [®]	913I-02532	1
Caelo	Unguentum Cordes [®]	9A163	1
Ichthyol	Unguentum Cordes [®]	11a042	1
Caelo	Unguentum Cordes [®]	913M-02532	1
Ichthyol / Noweda	Unguentum Cordes [®]	11A042	1
Ichthyol	Unguentum Cordes [®]	11A118	2
Ichthyol	Unguentum Cordes [®]	11A085	3
Ichthyol - Gesellschaft	Unguentum Cordes [®]	11A085	1
Ichthyol / Noweda	Unguentum Cordes [®]	11A146	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	11A177	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	11A171	1
Ichthyol/Anzag	Unguentum Cordes [®]	11A187	1
Kehr	Unguentum Cordes [®]	11A187	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	11A187	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	11A199	1
CIchthyol/K	Unguentum Cordes [®]	11A199	1
Ichthyol-Gesellschaft / K	Unguentum Cordes [®]	11A202	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	11A202	1
Ichthyol / Gehe	Unguentum Cordes [®]	11A202	1
Noweda	Unguentum Cordes [®]	11A227	1
Gehe	Unguentum Cordes [®]	11A223	2
Ichthyol	Unguentum Cordes [®]	11A233	1
IchthyolAnzag 10.05.2012	Unguentum Cordes [®]	11A234	1
Phönix Gotha	Unguentum Cordes [®]	11A234	1
Caelo	Unguentum Cordes [®]		1
Ichthyol/Otto	Unguentum Cordes [®]	12A020	1
Phönix	Unguentum Cordes [®]	12A049	1
Caelo	Unguentum Cordes [®]	12A022	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A054	1
Phoenix	Unguentum Cordes [®]	12A056	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A058	1
Ichtyol-Gesellschaft/Ph...	Unguentum Cordes [®]	12A058	1
Ichtyol 004	Unguentum Cordes [®]	12A059	1
Ichthyol/Gehe	Unguentum Cordes [®]	12A060	2
Ichthyol/Noweda	Unguentum Cordes [®]	12A064	1
Ichthyol-Gesellschaft / G	Unguentum Cordes [®]	12A064	1
Ichthyol	Unguentum Cordes [®]	12A078	1
Anzag	Unguentum Cordes [®]	12A077	1
Ichthyol/Gehe	Unguentum Cordes [®]	12A079	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	12A102	1
Cordes	Unguentum Cordes [®]	12A079	1
Ichthyol / Anzag	Unguentum Cordes [®]	12A079	1
Ichthyol / Gehe	Unguentum Cordes [®]	12A102	1
Ichthyol	Unguentum Cordes [®]	12A102	2
Ichthyol/Krieger	Unguentum Cordes [®]	12A102	1
Ichtyol	Unguentum Cordes [®]	12A106	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	12A115	2
Ichthyol / Gehe	Unguentum Cordes [®]	12A115	1
Ichthyol	Unguentum Cordes [®]	12A115	3
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	12A115	3
Ichthyol-Gesellschaft / G	Unguentum Cordes [®]	12A118	2
Ichthyol / Phoenix Leip...	Unguentum Cordes [®]	12A118	1
Ichthyol	Unguentum Cordes [®]	12A121	2
Ichthyol-Gesellschaft / G	Unguentum Cordes [®]	12A123	1
Cordes	Unguentum Cordes [®]	12A123	1
Noweda	Unguentum Cordes [®]	12A123	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	12A123	1
Bombastus	Unguentum Cordes [®]	12A123	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda/ Ichthyol	Unguentum Cordes [®]	12A124	1
Ichthyol/Anzag	Unguentum Cordes [®]	12A130	1
Ichthyol/Otto	Unguentum Cordes [®]	12A130	2
Sanacorp	Unguentum Cordes [®]	12A129	1
Anzag	Unguentum Cordes [®]	12A130	1
Sanova	Unguentum Cordes [®]	12A130	1
Phönix	Unguentum Cordes [®]	12A131	1
Ichthyol	Unguentum Cordes [®]	12A131	2
Phönix/Ichthyol	Unguentum Cordes [®]	12A139	1
Ichthyol/ Cordes	Unguentum Cordes [®]	12A139	1
Ichthyol	Unguentum Cordes [®]	12a145	1
Ichthtyol	Unguentum Cordes [®]	12a152	1
Caelo	Unguentum Cordes [®]	12A146	1
Corde	Unguentum Cordes [®]	12a152	1
Ichthyol	Unguentum Cordes [®]	12A152	2
IchthyolGesellschaft Co...	Unguentum Cordes [®]	12A152	1
ichthyol	Unguentum Cordes [®]	12A139	1
ICHTHYOL-GES/GEHE	Unguentum Cordes [®]	12A152	1
Maros / Phoenix	Unguentum Cordes [®]	12A152	1
ichthyol	Unguentum Cordes [®]	12A152	1
Ichthyol / Gehe	Unguentum Cordes [®]	12A152	2
Caelo	Unguentum Cordes [®]	12A152	1
Ichthyol /Noweda	Unguentum Cordes [®]	12A152	1
Fagron	Unguentum Cordes [®]	12A152	1
Ichthyol-Gesellschaft / E	Unguentum Cordes [®]	12A152	1
Ichthyol / Noweda	Unguentum Cordes [®]	12A152	2
Ichthyol	Unguentum Cordes [®]	12a152	1
Cordes	Unguentum Cordes [®]	12A155	1
Gehe	Unguentum Cordes [®]	12A155	1
Gehe	Unguentum Cordes [®]	12a155	2
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	12A155	1
Ichthyol/Gehe	Unguentum Cordes [®]	12A156	1
phönix/Ichthyol-Gesells...	Unguentum Cordes [®]	12A156	1
Caelo	Unguentum Cordes [®]	12A156	1
Ichtyol	Unguentum Cordes [®]	12A157	1
Sanacorp	Unguentum Cordes [®]	12A160	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	12A183	1
Ichthyol	Unguentum Cordes [®]	12A160	1
Caelo	Unguentum Cordes [®]	12A183	2
Ichthyol-Gesellschaft/V...	Unguentum Cordes [®]	12A183	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	12A183	1
Gehe	Unguentum Cordes [®]	12A183	1
IchthyolGesellschaft Co...	Unguentum Cordes [®]	12A183	1
Ichthyol / Sanacorp	Unguentum Cordes [®]	12A183	1
Ichthyol/ Noweda	Unguentum Cordes [®]	12A183	1
ichthyol/Sanacorp	Unguentum Cordes [®]	12A183	1
Noweda	Unguentum Cordes [®]	12A183	2
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	12A184	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	12A183	1
Gehe - Ichthyol	Unguentum Cordes [®]	12A183	1
Ichthyol / Phoenix	Unguentum Cordes [®]	12A183	1
Sanacorp	Unguentum Cordes [®]	12a184	1
Ichthyol/Krieger	Unguentum Cordes [®]	12A184	1
ICHTHYOL-Gesellschaft	Unguentum Cordes [®]	12A184	1
Ichthyol/Phönix	Unguentum Cordes [®]	12A184	1
Ichthyol-Ges	Unguentum Cordes [®]	12A184	1
Ichthyol / Gehe	Unguentum Cordes [®]	12A184	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	12A184	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A184	1
Ichthyol	Unguentum Cordes [®]	12A185	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol	Unguentum Cordes®	12A185	1
Ichthyol/AHD	Unguentum Cordes®	12A185	1
Caelo	Unguentum Cordes®	12A185	1
Ichthy	Unguentum Cordes®	12A185	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes®	12A185	1
Phoenix 26.11.12	Unguentum Cordes®	12A185	1
Stumpf/Ichthyol	Unguentum Cordes®	12A190	1
ichthyol	Unguentum Cordes®	12A190	1
Ichthyol	Unguentum Cordes®	12A190	1
Ichthyol/Gehe	Unguentum Cordes®	12A190	1
Gehe	Unguentum Cordes®	12A205	2
Ichthyol-Gesellschaft	Unguentum Cordes®	12A205	1
Ichthyol	Unguentum Cordes®	12a205	1
Ichthyol/ Cordes	Unguentum Cordes®	12A205	1
Ichthyol	Unguentum Cordes®	12A205	2
Caelo	Unguentum Cordes®	12A206	4
Ichthyol-Gesellschaft	Unguentum Cordes®	12A206	1
Ichthyol	Unguentum Cordes®	12A206	1
ichthyol	Unguentum Cordes®	12A206	1
Sanacorp / Ichthyol	Unguentum Cordes®	12A206	1
Phönix	Unguentum Cordes®	12A211	1
Ichthyol	Unguentum Cordes®	12A211	1
Ichthyol/Otto	Unguentum Cordes®	12AI30	1
Ichthyol/Sanacorp	Unguentum Cordes®	12A211	1
Ichthyol	Unguentum Cordes®	13092013D	1
Caelo	Unguentum Cordes®	131A043	1
Ichthyol	Unguentum Cordes®	13A005	1
Ichthyol	Unguentum Cordes®	13A003	1
Ichthyol/Gehe	Unguentum Cordes®	13A006	1
Ichthyol,Noweda	Unguentum Cordes®	13a009	2
Großhandel	Unguentum Cordes®	13A010	1
Ichthyol	Unguentum Cordes®	13A015	1
Ichthyol	Unguentum Cordes®	13A006	1
Cordes	Unguentum Cordes®	13a025	1
IchthyolSanacorp	Unguentum Cordes®	13A025	1
Ichthyol	Unguentum Cordes®	13A010	2
Ichthyol	Unguentum Cordes®	13A029	2
Ichthyol/Phönix	Unguentum Cordes®	13A029	1
Anzag/ Ichtyol	Unguentum Cordes®	13A029	1
Ichthyol/Phönix	Unguentum Cordes®	13A030	1
Ichthyol	Unguentum Cordes®	13A030	2
Ichthyol	Unguentum Cordes®	13A031	2
Ichthyol/Phönix	Unguentum Cordes®	13A031	1
Gehe	Unguentum Cordes®	13A031	3
Alliance Health	Unguentum Cordes®	13A032	1
ICHTHYOL/11	Unguentum Cordes®	13A031	1
ichthyol	Unguentum Cordes®	13A032	1
Ichthyol	Unguentum Cordes®	13A032	2
Phönix Gotha	Unguentum Cordes®	13A032	1
Ichthyol/Phönix	Unguentum Cordes®	13A033	1
Phönix	Unguentum Cordes®	13A034	1
Ichthyol/phönix	Unguentum Cordes®	13A034	1
Ichthyol	Unguentum Cordes®	13A034	2
Ichthyol-Gesellschaft, ...	Unguentum Cordes®	13A034	1
Ichthyol-Gesellschaft	Unguentum Cordes®	13A034	1
Ichthyol/GEHE	Unguentum Cordes®	13A034	1
Ichthyol/ebert	Unguentum Cordes®	13A034	1
Ichthyol/ Noweda	Unguentum Cordes®	13A034	1
	Unguentum Cordes®	13A034	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	13A034	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/Noweda	Unguentum Cordes [®]	13A042	1
Ichthyol	Unguentum Cordes [®]	13A042	2
Ichthyol / Noweda	Unguentum Cordes [®]	13A042	1
Gehe/Ichthyol	Unguentum Cordes [®]	13A042	1
Bombastus	Unguentum Cordes [®]	13A042	1
Ichthyol / Gehe	Unguentum Cordes [®]	13A042	1
Ichthyol-Gesellschaft / G	Unguentum Cordes [®]	13A042	1
Ichthyol/Noweda	Unguentum Cordes [®]	13A042	1
Ichthyol	Unguentum Cordes [®]	13A043	2
Ichthyol / Noweda	Unguentum Cordes [®]	13A043	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	13A043	1
Ichthyol	Unguentum Cordes [®]	13A043	2
Caelo	Unguentum Cordes [®]	13a043	1
Caelo	Unguentum Cordes [®]	13A043	1
Ichthyol/ Noweda	Unguentum Cordes [®]	13A043	1
Ichthyol/Alliance healt...	Unguentum Cordes [®]	13A047	1
Ichthyol	Unguentum Cordes [®]	13A047	2
Ichthyol/Noweda	Unguentum Cordes [®]	13A043	1
Gehe	Unguentum Cordes [®]	13A047	1
ichthyol/ANZAG	Unguentum Cordes [®]	13A047	1
cordes/phoenix	Unguentum Cordes [®]	13A047	1
Caelo	Unguentum Cordes [®]	13A049	1
Ichthyol	Unguentum Cordes [®]	13A048	1
Ichthyol/Phoenix	Unguentum Cordes [®]	13A049	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	13A048	1
Jenne	Unguentum Cordes [®]	13A049	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	13A050	1
Ichthyol/Anzag	Unguentum Cordes [®]	13A049/13A050	1
Ichthyol-Gesellschaft, ...	Unguentum Cordes [®]	13A050	1
Cordes/Phönix	Unguentum Cordes [®]	13A050	1
Ichthyol	Unguentum Cordes [®]	13A050	4
Ichthyol/Gehe	Unguentum Cordes [®]	13A051	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	13A051	1
Stumpf/Ichthyol	Unguentum Cordes [®]	13A051	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	13A051	1
phönix	Unguentum Cordes [®]	13A051	1
Ichthyol / Noweda	Unguentum Cordes [®]	13A051	1
Ichthyol	Unguentum Cordes [®]	13A052	1
ichthyol gesellschaft	Unguentum Cordes [®]	13A051	1
Cordes	Unguentum Cordes [®]	13A052	1
Ebert/Ichthyol	Unguentum Cordes [®]	13A052	1
Hedinger	Unguentum Cordes [®]	13A053	1
Ichthyol/Noweda	Unguentum Cordes [®]	13A053	1
Ichthyol	Unguentum Cordes [®]	13A053	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	13A053	1
Gehe	Unguentum Cordes [®]	13A054	1
Ichthyol-Gesellschaft/P...	Unguentum Cordes [®]	13A053	1
Ichthyol / Gehe	Unguentum Cordes [®]	13A054	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	13A054	1
Ichthyol/Phoenix	Unguentum Cordes [®]	13A054	1
Ichthyol/Otto	Unguentum Cordes [®]	13A055	1
Caelo	Unguentum Cordes [®]	13A055	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	13A055	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	13A055	1
Ichthyol	Unguentum Cordes [®]	13A055	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	13A055	1
Ichthyol/Noweda	Unguentum Cordes [®]	13A055	1
phönix/Ichthyol-Gesells...	Unguentum Cordes [®]	13A055	1
Ichthyol/Anzag	Unguentum Cordes [®]	13A067	1
Ichthyol	Unguentum Cordes [®]	13A067	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
ichthyol gesellschaft	Unguentum Cordes [®]	13A055	1
Ichthyol	Unguentum Cordes [®]	13A068	2
Ichthyol Gesellschaft/G...	Unguentum Cordes [®]	13A089	1
Gehe/Ichthyol	Unguentum Cordes [®]	13A068	1
Ichthyol	Unguentum Cordes [®]	13A089	2
Ichthyol	Unguentum Cordes [®]	13A067	1
Sanacorp	Unguentum Cordes [®]	13A068	1
Euro OTC	Unguentum Cordes [®]	13A068	1
Ichthyol /Noweda	Unguentum Cordes [®]	13A090	2
Ichthyol	Unguentum Cordes [®]	13A090	1
Cordes	Unguentum Cordes [®]	13A106	1
Ichthyol/Noweda	Unguentum Cordes [®]	13A106	1
Noweda/Ichthyol	Unguentum Cordes [®]	13A106	1
Gehe	Unguentum Cordes [®]	13A106	1
Ichthyol / Gehe	Unguentum Cordes [®]	13A106	1
Ichthyol/Gehe	Unguentum Cordes [®]	13A106	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	13A106	1
Ichthyol	Unguentum Cordes [®]	13A106	2
Ichthyol	Unguentum Cordes [®]	13A107	1
Ichthyol	Unguentum Cordes [®]	13a107	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	13A108	1
Ichthyol/Gehe	Unguentum Cordes [®]	13A109	1
Ichthiol	Unguentum Cordes [®]	13A110	1
Ichthyol	Unguentum Cordes [®]	13A111	1
ichthyol/Gehe	Unguentum Cordes [®]	13A109	1
Ichtyol	Unguentum Cordes [®]	13A111	1
Ichthyol/Phoenix	Unguentum Cordes [®]	13A135	1
Ichthyol	Unguentum Cordes [®]	13A112	2
Sanacorp	Unguentum Cordes [®]	13A136	1
Ichtyol/ Phönix	Unguentum Cordes [®]	13A135	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	13A135	1
Gehe	Unguentum Cordes [®]	13A136	2
Ichthyol/VDL	Unguentum Cordes [®]	13A136	1
ICHTHYOL	Unguentum Cordes [®]	13A135	1
Ichthyol/Otto	Unguentum Cordes [®]	13A136	1
Ichthyol	Unguentum Cordes [®]	13A136	1
Ichtyol	Unguentum Cordes [®]	13A137	1
Gehe - Ichthyol	Unguentum Cordes [®]	13A136	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	13A137	1
Ichthyol GmbH	Unguentum Cordes [®]	13A137	1
Caelo	Unguentum Cordes [®]	13A138	2
Ichthyol / Gehe	Unguentum Cordes [®]	13A138	1
Gehe	Unguentum Cordes [®]	13A141	1
Gehe	Unguentum Cordes [®]	13A153	1
Stumpf/Ichtyol	Unguentum Cordes [®]	13A142	1
Stumpf/Ichtyol	Unguentum Cordes [®]	13A166	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	13A166	1
Ichthtyol	Unguentum Cordes [®]	13A178	1
Ichthyol/Krieger/23.05...	Unguentum Cordes [®]	13A209	1
noweda/Ichtyol	Unguentum Cordes [®]	13A043	1
Ichtyol,Noweda	Unguentum Cordes [®]	13A047	1
Ichthyol/Phönix	Unguentum Cordes [®]	13A054	1
Caelo	Unguentum Cordes [®]	13A047	1
ichthyol	Unguentum Cordes [®]	1413M-02532	1
Ichthyol	Unguentum Cordes [®]	14A005	1
Ichtyol	Unguentum Cordes [®]	13A054	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	14A055	2
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	14A056	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	14A060	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	14A060	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	14A113	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	14A116	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	14A130	1
ichthyol gesellschaft	Unguentum Cordes [®]	14A143	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	14A168	1
Caelo	Unguentum Cordes [®]	14A143	1
ichthyol gesellschaft	Unguentum Cordes [®]	15A022	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A038	1
ichthyol gesellschaft	Unguentum Cordes [®]	15A041	1
Cordes	Unguentum Cordes [®]	15A060	1
Ichthyol	Unguentum Cordes [®]	15A075	1
ichthyol gesellschaft	Unguentum Cordes [®]	15A076	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A076	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A079	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A080	1
Ichthyol	Unguentum Cordes [®]	15A123	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A139	1
ichtyol	Unguentum Cordes [®]	15a195	1
ichthyol gesellschaft	Unguentum Cordes [®]	15A159	2
ichthyol gesellschaft	Unguentum Cordes [®]	16A011	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A073	1
Caelo	Unguentum Cordes [®]	16A074	1
Sanacorp	Unguentum Cordes [®]	2/AKD2	1
ichthyol	Unguentum Cordes [®]	20201112-1	1
Ichthyol	Unguentum Cordes [®]	24022014B	1
Caelo	Unguentum Cordes [®]	3112A-02532	1
Ichtyol/Phönix	Unguentum Cordes [®]	4113Q-02532	1
Caelo	Unguentum Cordes [®]	Ichtyol	1
Ichthyol	Unguentum Cordes [®]	11A202	2
Lamotte/Noweda	Paraffin, liquid	1000002995	1
Gehe/Hedinger	Paraffin, liquid	1000002995	3
Fagron	Paraffin, liquid	16E30-B05-327577	4
Caelo	Paraffin, liquid	17062005	2
Caelo	Paraffin, liquid	17062007	4
Caelo	Paraffin, liquid	1701033-01	1
Caelo	Paraffin, liquid	17321401	1
Sanacorp WE: 20.11.17 1...	Paraffin, liquid	17062002	1
Phönix	Paraffin, liquid	17062007	1
Fagron	Paraffin, liquid	18B19-B06	1
Caelo	Paraffin, liquid	17047603	1
Apomix	Eye ointment, emulsifying DAC	2B020117	1
Gehe/Apomix	Eye ointment, emulsifying DAC	1C080114	1
Apomix	Eye ointment, emulsifying DAC	1D040116	1
Apomix /Noweda	Eye ointment, emulsifying DAC	2D020117	1
Apomix/Kehr	Eye ointment, emulsifying DAC	2D040116	1
Apomix/Phönix	Eye ointment, emulsifying DAC	3C060217	1
Apomix / Kehr	Eye ointment, emulsifying DAC	48200114	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	4B200114	1
Noweda/apomix	Eye ointment, emulsifying DAC	4B200114	1
Caelo	Eye ointment, emulsifying DAC	5D110416	1
Apomix	Eye ointment, emulsifying DAC	5A030314	1
Apomix / Noweda	Eye ointment, emulsifying DAC	6A090516	1
AMH Niemannn GmbH	Eye ointment, emulsifying DAC	5B110416	1
Apomix/Kehr	Eye ointment, emulsifying DAC	7A160315	1
Allianc. Health. /Apomix	Eye ointment, emulsifying DAC	5D110416	1
Apomix/Kehr	Eye ointment, emulsifying DAC	7D180417	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	6D150317	1
apomix/Holdermann	Eye ointment, emulsifying DAC	7A180417	1
	Eye ointment, emulsifying DAC	68150317	1
AHD / apomix	Eye ointment, emulsifying DAC	7B010616	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix Noweda	Eye ointment, emulsifying DAC	8B150517	1
Apomix/Noweda	Eye ointment, emulsifying DAC	8D020414	1
apomix/Sanacorp	Eye ointment, emulsifying DAC	9A190514	1
Apomix / Noweda	Eye ointment, emulsifying DAC	9A080415	1
Apomix/Anzag	Eye ointment, emulsifying DAC	9A080415	1
Apomix	Eye ointment, emulsifying DAC	9C050916	1
Apomix/Noweda	Eye ointment, emulsifying DAC	10A100717	1
Alliance Healthcare/Apo...	Eye ointment, emulsifying DAC	9C190514	1
Allianc. Healthc. /Apomix	Eye ointment, emulsifying DAC	10A010615	2
Allianc. Health. / Apomix	Eye ointment, emulsifying DAC	10C041016	2
Caelo	Emulsifying hydrophobic base...	13340814	1
Caelo	Emulsifying hydrophobic base...	13295801	3
Apomix/Phönix	Eye ointment, emulsifying DAC	11A011116	1
Caelo	Emulsifying hydrophobic base...	13361004	9
Caelo	Emulsifying hydrophobic base...	14087501	7
Caelo	Emulsifying hydrophobic base...	14161501	11
Caelo	Emulsifying hydrophobic base...	14276301	5
Caelo	Emulsifying hydrophobic base...	14276401	12
Alliance Healthcare/Apo...	Eye ointment, emulsifying DAC	14A011013	1
Caesar & Loretz GmbH	Emulsifying hydrophobic base...	14276401	2
Apomix / Kehr	Eye ointment, emulsifying DAC	14A011013	1
Apomix	Eye ointment, emulsifying DAC	14C170914	1
Noweda	Emulsifying hydrophobic base...	14276401	2
Phönix Gotha	Emulsifying hydrophobic base...	15064101	1
Caelo	Emulsifying hydrophobic base...	15064101	12
Caelo	Emulsifying hydrophobic base...	15205002	13
Caelo	Emulsifying hydrophobic base...	15269201	13
Caelo	Emulsifying hydrophobic base...	15331102	2
Caelo	Emulsifying hydrophobic base...	15331101	2
Caelo	Emulsifying hydrophobic base...	16054201	8
Caelo	Emulsifying hydrophobic base...	16054302	7
Euro OTC	Emulsifying hydrophobic base...	16172201	3
Caelo	Emulsifying hydrophobic base...	16172201	3
Caelo	Emulsifying hydrophobic base...	16172301	6
Caelo	Emulsifying hydrophobic base...	16172401	4
Caelo	Emulsifying hydrophobic base...	16172402	1
Caelo	Emulsifying hydrophobic base...	16225902	1
Caelo	Emulsifying hydrophobic base...	17032201	4
Gehe	Emulsifying hydrophobic base...	17032303	1
Caelo	Emulsifying hydrophobic base...	17032303	4
Noweda	Emulsifying hydrophobic base...	17113703	1
Noweda	Emulsifying hydrophobic base...	17032303	1
Caelo	Emulsifying hydrophobic base...	17113703	3
apomix/Noweda	Eye ointment, emulsifying DAC	18C021115	1
Alliance HealthcareApomix	Eye ointment, emulsifying DAC	18B021214	1
AHD	Eye ointment, emulsifying DAC	18A021115	1
Allianc. Health. /Apomix	Eye ointment, emulsifying DAC	18A021115	1
Noweda/apomix	Eye ointment, emulsifying DAC	18C021214	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	18C021214	1
PKH/Gehe	Eye ointment, emulsifying DAC	3B180213	2
Apomix/Kehr	Eye ointment, emulsifying DAC	20D011215	1
Caelo	Eye ointment, emulsifying DAC	4A270212	1
Apomix/Noweda	Eye ointment, emulsifying DAC	7D210513	1
Apomix	Eye ointment, emulsifying DAC	7C040612	1
PKH Halle / Noweda	Eye ointment, emulsifying DAC	9C250612	1
Apomix / Sanacorp	Eye ointment, emulsifying DAC	12C011012	1
PKH/Gehe	Eye ointment, emulsifying DAC	10C170613	2
apomix/ Sanacorp	Eye ointment, emulsifying DAC	12C300813	1
Apomix / AHD	Eye ointment, emulsifying DAC	11B170713	1
Noweda/ PKH Halle	Eye ointment, emulsifying DAC	16A061212	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix/Anzag	Eye ointment, emulsifying DAC	16B061212	1
Caelo	Emulsifying hydrophobic base...	12182301	1
Gehe 07.03.13	Emulsifying hydrophobic base...	271904	2
Caelo	Emulsifying hydrophobic base...	12229003	4
Caelo	Emulsifying hydrophobic base...	12343501	6
Caelo	Emulsifying hydrophobic base...	13066802	10
gehe 13.07.2013	Emulsifying hydrophobic base...	13066802	2
Phönix 25.09.13	Emulsifying hydrophobic base...	13165201	1
Phoenix 05.11.2013	Emulsifying hydrophobic base...	13165205	1
Caelo	Emulsifying hydrophobic base...	13177601	7
Caelo	Emulsifying hydrophobic base...	13295701	8
Caelo	Emulsifying hydrophobic base...	20081302	1
Bombastus	Emulsifying hydrophobic base...	279074	1
Bombastus	Emulsifying hydrophobic base...	22011402	1
Bombastus	Wool alcohol ointment DAB	915M-01370	1
Beiersdorf AG	Wool alcohol ointment DAB	11052015B	1
Caelo	Wool alcohol ointment DAB	14282901	5
Caelo	Wool alcohol ointment DAB	15035503	3
Caelo	Wool alcohol ointment DAB	15228202	16
Caelo	Wool alcohol ointment DAB	16007201	6
Caelo	Wool alcohol ointment DAB	15361801	11
Caelo	Wool alcohol ointment DAB	15225601	11
Caesar & Loretz/Noweda	Wool alcohol ointment DAB	16066601	1
Caelo	Wool alcohol ointment DAB	16066601	6
EuroOTV	Wool alcohol ointment DAB	1608040-01	1
Caelo	Wool alcohol ointment DAB	160905EA	1
Sanacorp	Wool alcohol ointment DAB	28/AKD12	1
Bombastus	Wool alcohol ointment DAB	293093	6
Bombastus	Wool alcohol ointment DAB	298579	3
Sanacorp	Wool alcohol ointment DAB		1
Infectopharm/ Sanacorp	Wool alcohol ointment DAB	50334610	1
Caelo	Wool alcohol ointment DAB	50815310	1
Beiersdorf AG/Sanacorp	Wool alcohol ointment DAB	50444710	1
Beiersdorf	Wool alcohol ointment DAB	50815310	1
Beiersdorf / Phoenix	Wool alcohol ointment DAB	50845210	1
Caelo	Wool alcohol ointment DAB	50845210	1
Noweda	Wool alcohol ointment DAB	50845210	2
Beiersdorf/Holdermann	Wool alcohol ointment DAB	50845210	1
Beiersdorf AG/NOWEDA	Wool alcohol ointment DAB	50845210	1
Beiersdorf	Wool alcohol ointment DAB	50855310	1
Beiersdorf/N	Wool alcohol ointment DAB	50925410	1
Beiersdorf	Wool alcohol ointment DAB	50925410	6
Beiersdorf / Noweda	Wool alcohol ointment DAB	50925410	1
Beiersdorf	Wool alcohol ointment DAB	50935410	3
Beiersdorf / Phoenix	Wool alcohol ointment DAB	50935410	1
Beiersdorf/Beiersdorf	Wool alcohol ointment DAB	50955510	1
Schärer & Schlöpfer AG/...	Wool alcohol ointment DAB	50955510	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	51015610	1
Beiersdorf AG	Wool alcohol ointment DAB	51025610	1
Beiersdorf /Sanacorp	Wool alcohol ointment DAB	51025610	1
Beiersdorf	Wool alcohol ointment DAB	51015610	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	51115710	2
Beiersdorf	Wool alcohol ointment DAB	51125710	2
Beiersdorf	Wool alcohol ointment DAB	51115710	1
Beiersdorf	Wool alcohol ointment DAB	51946010	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	51946010	3
Beiersdorf direkt	Wool alcohol ointment DAB	52056110	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	52656810	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	53437010	2
Caelo	Wool alcohol ointment DAB	53437010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Beiersdorf/N	Wool alcohol ointment DAB	53437010	1
Beiersdorf	Wool alcohol ointment DAB	53437010	2
Beiersdorf/Gehe	Wool alcohol ointment DAB	53447010	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	53437010	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	53447010	1
Caelo	Wool alcohol ointment DAB	53457110	1
Beiersdorf/N	Wool alcohol ointment DAB	53457110	1
Beiersdorf / Sanacorp	Wool alcohol ointment DAB	53457110	1
Beiersdorf/Sana	Wool alcohol ointment DAB	54137710	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	54147810	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	54147810	2
Beiersdorf /Alliance He...	Wool alcohol ointment DAB	54127710	1
Beiersdorf	Wool alcohol ointment DAB	54127710	2
Beiersdorf/ Phoenix	Wool alcohol ointment DAB	54127710	1
Beiersdorf/N	Wool alcohol ointment DAB	54147810	1
Beiersdorf	Wool alcohol ointment DAB	54147810	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	54157810	1
Beiersdorf	Wool alcohol ointment DAB	54157810	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	54157810	1
Caelo	Wool alcohol ointment DAB	54157810	1
Beiersdorf	Wool alcohol ointment DAB	60419110	3
Beiersdorf/E	Wool alcohol ointment DAB	60419110	1
Gehe	Wool alcohol ointment DAB	60419110	1
Noweda/Beiersdorf	Wool alcohol ointment DAB	60539210	1
Beiersdorf	Wool alcohol ointment DAB	60539210	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	60839310	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	61079710	1
Beiersdorf/N	Wool alcohol ointment DAB	60839310	1
Beiersdorf AG	Wool alcohol ointment DAB	60919510	2
Beiersdorf	Wool alcohol ointment DAB	61079710	1
Beiersdorf / Kehr	Wool alcohol ointment DAB	61149910	1
Beiersdorf AG	Wool alcohol ointment DAB	61149910	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	61841610	1
Beiersdorf AG	Wool alcohol ointment DAB	61841610	1
Beiersdorf / Phoenix	Wool alcohol ointment DAB	62031710	2
Beiersdorf	Wool alcohol ointment DAB	62031710	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	62031710	1
Ebert 18.10.16	Wool alcohol ointment DAB	62031710	1
Allianc. Health. /Beier...	Wool alcohol ointment DAB	62031710	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	62041710	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	62922210	1
Caelo	Wool alcohol ointment DAB	62922210	1
Beiersdorf/AHC	Wool alcohol ointment DAB	62932210	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	62932210	1
Beiersdorf	Wool alcohol ointment DAB	62932210	1
Beiersdorf/Anzag	Wool alcohol ointment DAB	63052310	1
Beiersdorf AG	Wool alcohol ointment DAB	63052310	1
Caelo	Wool alcohol ointment DAB	63052310	2
Beiersdorf/Phönix	Wool alcohol ointment DAB	64312910	2
Beiersdorf AG	Wool alcohol ointment DAB	64312910	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	64332810	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	64342810	1
Gehe/Beiersdorf	Wool alcohol ointment DAB	64342810	1
Beiersdorf	Wool alcohol ointment DAB	64342810	2
Caelo	Wool alcohol ointment DAB	64342810	2
Beiersdorf AG	Wool alcohol ointment DAB	64713510	1
Beiersdorf	Wool alcohol ointment DAB	64713510	3
beiersdorf	Wool alcohol ointment DAB	64833610	1
Beiersdorf/ Noweda	Wool alcohol ointment DAB	64833610	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	70314010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Beiersdorf/Anzag	Wool alcohol ointment DAB	71114510	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	71114510	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	71114510	1
Beiersdorf / Phoenix	Wool alcohol ointment DAB	71114510	1
Caelo	Wool alcohol ointment DAB	71124510	1
Noweda/Beiersdorf	Wool alcohol ointment DAB	71214810	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	72915810	1
Beiersdorf AG	Wool alcohol ointment DAB	71274810	1
Caelo	Wool alcohol ointment DAB	74017010	1
Beiersdorf/Kehr	Wool alcohol ointment DAB	72915810	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	74747510	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	74747510	1
Noweda/Beiersdorf	Wool alcohol ointment DAB	74757510	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	81819410	1
Beiersdorf AG	Wool alcohol ointment DAB	10112014A	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	7477510	1
Bombastus	Wool alcohol ointment DAB	1211018-01	1
Beiersdorf AG	Wool alcohol ointment DAB	13042015D	1
Caelo	Wool alcohol ointment DAB	131255101	1
Caelo	Wool alcohol ointment DAB	13255001	10
Caelo	Wool alcohol ointment DAB	13242602	15
Caelo	Wool alcohol ointment DAB	13242601	13
Sanacorp	Wool alcohol ointment DAB	13255001	1
Caelo	Wool alcohol ointment DAB	13255101	24
Caelo	Wool alcohol ointment DAB	13255201	10
Caelo	Wool alcohol ointment DAB	14067002	6
Caelo	Wool alcohol ointment DAB	14153101	12
Caelo	Wool alcohol ointment DAB	14067201	13
Caelo	Wool alcohol ointment DAB	14153602	3
Caelo	Wool alcohol ointment DAB	14153603	8
Caelo	Wool alcohol ointment DAB	14282801	8
Phönix Gotha	Wool alcohol ointment DAB	14282801	1
Caelo	Wool alcohol ointment DAB	14283101	3
Caelo	Wool alcohol ointment DAB	15034801	7
Caelo	Wool alcohol ointment DAB	15035301	9
Caelo	Wool alcohol ointment DAB	15034901	5
Bombastus	Wool alcohol ointment DAB	289510	3
Bombastus	Wool alcohol ointment DAB	279927	3
Bombastus	Wool alcohol ointment DAB	287012	1
Beiersdorf AG	Wool alcohol ointment DAB	33753710	1
Beiersdorf	Wool alcohol ointment DAB	34554710	5
Beiersdorf/Phönix	Wool alcohol ointment DAB	34554710	1
Anzag/Beiersdorf	Wool alcohol ointment DAB	40930610	1
Beiersdorf	Wool alcohol ointment DAB	40726610	1
Stumpf/Beiersdorf Ag	Wool alcohol ointment DAB	40956810	1
Großhandel	Wool alcohol ointment DAB	40956810	1
Beiersdorf AG	Wool alcohol ointment DAB	40956810	1
Beiersdorf / Sanacorp	Wool alcohol ointment DAB	41026910	1
Beiersdorf AG / Phönix	Wool alcohol ointment DAB	40956810	1
Beiersdorf	Wool alcohol ointment DAB	40956810	2
Beiersdorf	Wool alcohol ointment DAB	41026910	2
Gehe	Wool alcohol ointment DAB	41026910	1
Caelo	Wool alcohol ointment DAB	41026910	1
Beiersdorf/Kehr	Wool alcohol ointment DAB	4104701008/16	1
Beiersdorf	Wool alcohol ointment DAB	41047010	2
Beiersdorf	Wool alcohol ointment DAB	41839810	2
Stumpf/Beiersdorf Ag	Wool alcohol ointment DAB	41910710	1
Beiersdorf/Otto	Wool alcohol ointment DAB	41047010	1
Beiersdorf AG	Wool alcohol ointment DAB	41829810	1
Sanacorp	Wool alcohol ointment DAB	41047010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Beiersdorf/Noweda	Wool alcohol ointment DAB	41910710	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	41910710	1
Noweda / Beiersdorf AG	Wool alcohol ointment DAB	41930610	1
Beiersdorf AG	Wool alcohol ointment DAB	41930610	2
Beiersdorf/Gehe	Wool alcohol ointment DAB	41930610	1
Beiersdorf	Wool alcohol ointment DAB	41940610	1
Gehe/ Beiersdorf AG	Wool alcohol ointment DAB	41940610	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	41930610	1
Beiersdorf AG / Phönix	Wool alcohol ointment DAB	41940610	1
Beiersdorf/ Gehe	Wool alcohol ointment DAB	41930610	1
Beiersdorf/ Gehe	Wool alcohol ointment DAB	41950710	1
Beiersdorf SA	Wool alcohol ointment DAB	41950710	1
Beiersdorf	Wool alcohol ointment DAB	42531910	1
Sanacorp	Wool alcohol ointment DAB	43112910	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	43142910	2
BDF	Wool alcohol ointment DAB	43142910	1
Beiersdorf SA	Wool alcohol ointment DAB	44824310	1
Beiersdorf/ Noweda	Wool alcohol ointment DAB	43142910	1
Gehe	Wool alcohol ointment DAB	43142910	1
Gehe	Wool alcohol ointment DAB	45054410	1
Beiersdorf	Wool alcohol ointment DAB	45054410	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	45054410	2
Noweda / Beiersdorf AG	Wool alcohol ointment DAB	45054410	1
Caelo	Wool alcohol ointment DAB	50324610	1
Caelo	Wool alcohol ointment DAB	46/AKD2	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	50324610	1
Noweda/Beiersdorf	Wool alcohol ointment DAB	50324610	1
Beiersdorf	Wool alcohol ointment DAB	50324610	2
Caelo	Wool alcohol ointment DAB	50414810	1
Beiersdorf über Sanacorp	Wool alcohol ointment DAB	50444710	1
Beiersdorf AG	Wool alcohol ointment DAB	50444710	2
Gehe/Beiersdorf	Wool alcohol ointment DAB	50414810	1
Beiersdorf	Wool alcohol ointment DAB	50444710	4
Alliance Healthcare/Bei. . .	Wool alcohol ointment DAB	50454710	1
Noweda / Beiersdorf AG	Wool alcohol ointment DAB	50815310	1
Beiersdorf AG/Noweda	Wool alcohol ointment DAB	50815310	1
Beiersdorf AG	Wool alcohol ointment DAB	51946010	1
Beiersdorf AG	Wool alcohol ointment DAB	54738610	1
Beiersdorf	Wool alcohol ointment DAB	60919510	2
Beiersdorf	Wool alcohol ointment DAB	62922210	1
Baiersdorf	Wool alcohol ointment DAB	4833010	1
Beiersdorf AG	Wool alcohol ointment DAB	6022014A	2
Caelo	Wool alcohol ointment DAB	9365602	1
Caelo	Wool alcohol ointment DAB	12152001	5
Sanacorp	Wool alcohol ointment DAB	12151801	1
Bombastus	Wool alcohol ointment DAB	1109210014	1
Caelo	Wool alcohol ointment DAB	1212M-01370	1
Caelo	Wool alcohol ointment DAB		10
Fagron	Wool alcohol ointment DAB	1234	2
Caelo	Wool alcohol ointment DAB	13067102	14
Caelo	Wool alcohol ointment DAB	12340001	4
Caelo	Wool alcohol ointment DAB	13126102	1
Caelo	Wool alcohol ointment DAB	13114102	8
Phönix	Wool alcohol ointment DAB	13114103	1
Caelo	Wool alcohol ointment DAB	13114001	7
Caelo	Wool alcohol ointment DAB	13242901	5
Gehe	Wool alcohol ointment DAB	13357910	1
Beiersdorf (DAB 2010)	Wool alcohol ointment DAB	13658410	1
Beiersdorf	Wool alcohol ointment DAB	14021301	1
Caelo	Wool alcohol ointment DAB	14024701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Beiersdorf	Wool alcohol ointment DAB	14939610	1
Beiersdorf	Wool alcohol ointment DAB	20440810	1
Gehe	Wool alcohol ointment DAB	21241410	2
Beiersdorf	Wool alcohol ointment DAB	21241410	1
Beiersdorf	Wool alcohol ointment DAB	21431810	1
Beiersdorf (DAB 2010)	Wool alcohol ointment DAB	22412510	1
Beiersdorf	Wool alcohol ointment DAB	22142410	1
Gehe	Wool alcohol ointment DAB	22412510	2
Anzag	Wool alcohol ointment DAB	22412510	2
BDF Phönix	Wool alcohol ointment DAB	22422510	1
Phoenix	Wool alcohol ointment DAB	22532710	1
Noweda	Wool alcohol ointment DAB	22532710	1
Beiersdorf AG / Noweda	Wool alcohol ointment DAB	22532710	1
Phönix/Beiersdorf	Wool alcohol ointment DAB	23023210	1
Caelo	Wool alcohol ointment DAB	2313E-01370	1
Caelo	Wool alcohol ointment DAB	23023210	2
Stumpf	Wool alcohol ointment DAB	23223410	1
Beiersdorf AG	Wool alcohol ointment DAB	23223410	1
BDF	Wool alcohol ointment DAB	2322410	1
Beiersdorf	Wool alcohol ointment DAB	23223410	4
Anzag/	Wool alcohol ointment DAB	23223410	1
Beiersdorf	Wool alcohol ointment DAB	23753710	1
Gehe	Wool alcohol ointment DAB	23753710	1
Beiersdorf	Wool alcohol ointment DAB	24113810	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	24123810	1
Beiersdorf	Wool alcohol ointment DAB	24233910	2
Beiersdorf / Noweda	Wool alcohol ointment DAB	24233910	1
Phönix/Beiersdorf	Wool alcohol ointment DAB	24243910	1
Sanacorp	Wool alcohol ointment DAB	24254010	1
Anzag	Wool alcohol ointment DAB	24254010	1
Beiersdorf	Wool alcohol ointment DAB	24254010	1
Beiersdorf AG	Wool alcohol ointment DAB	24635410	1
Beiersdorf AG	Wool alcohol ointment DAB	24254010	1
Beiersdorf	Wool alcohol ointment DAB	24635410	4
Noweda	Wool alcohol ointment DAB	24635410	1
Bombastus	Wool alcohol ointment DAB	281690	5
Gehe	Wool alcohol ointment DAB	30711310	1
beiersdorf/Gehe	Wool alcohol ointment DAB	30721310	1
Beiersdorf Manufacturing	Wool alcohol ointment DAB	30731210	1
Beiersdorf S.A.	Wool alcohol ointment DAB	30731210	1
Caelo	Wool alcohol ointment DAB	30731210	1
Beiersdorf AG	Wool alcohol ointment DAB	30731210	2
Beiersdorf AG	Wool alcohol ointment DAB	30821610	1
Gehe	Wool alcohol ointment DAB	30831410	2
Beiersdorf	Wool alcohol ointment DAB	30831410	1
Beiersdorf AG/NOWEDA	Wool alcohol ointment DAB	30831410	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	30831410	1
Beiersdorf	Wool alcohol ointment DAB	30841410	4
Caelo	Wool alcohol ointment DAB	30841410	1
Stumpf	Wool alcohol ointment DAB	30851510	1
Beiersdorf S.A.	Wool alcohol ointment DAB	30851510	1
Sanacorp	Wool alcohol ointment DAB	30851510	1
Caelo	Wool alcohol ointment DAB	31121710	2
Beiersdorf AG	Wool alcohol ointment DAB	31231810	2
Gehe	Wool alcohol ointment DAB	30851510	1
alliance Healthcare/Bei...	Wool alcohol ointment DAB	30851510	1
Gehe	Wool alcohol ointment DAB	31231810	1
Phönix	Wool alcohol ointment DAB	31322010	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	31231810	1
Beiersdorf/Otto	Wool alcohol ointment DAB	31251910	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Beiersdorf/Noweda	Wool alcohol ointment DAB	31322010	2
Beiersdorf S.A.	Wool alcohol ointment DAB	32352710	1
Stumpf/Beiersdorf Ag	Wool alcohol ointment DAB	33323210	1
Beiersdorf/ Noweda	Wool alcohol ointment DAB	33323210	2
Beiersdorf	Wool alcohol ointment DAB	33323210	2
Beiersdorf AG/Noweda	Wool alcohol ointment DAB	33323210	2
Gehe	Wool alcohol ointment DAB	33323210	1
Beiersdorf/Krieger/06.1...	Wool alcohol ointment DAB	33323210	1
Beiersdorf/Gumo	Wool alcohol ointment DAB	33323210	1
Beiersdorf	Wool alcohol ointment DAB	33343310	3
BDF	Wool alcohol ointment DAB	33343310	1
Stumpf/Beiersdorf Ag	Wool alcohol ointment DAB	33343310	1
Beiersdorf AG	Wool alcohol ointment DAB	33343310	1
Beiersdorf	Wool alcohol ointment DAB	33723610	2
Beiersdorf/Phönix	Wool alcohol ointment DAB	33343310	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	33723610	1
Beiersdorf/V.d.L.	Wool alcohol ointment DAB	33723610	1
Gehe	Wool alcohol ointment DAB	33753710	3
Beiersdorf AG / Noweda	Wool alcohol ointment DAB	33753710	1
Caelo	Wool alcohol ointment DAB	33753710	1
Sanacorp/ Beiersdorf	Wool alcohol ointment DAB	33723610	1
Beiersdorf/Otto	Wool alcohol ointment DAB	33753710	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	33753710	2
Beiersdorf	Wool alcohol ointment DAB	34224110	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	34254110	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	34414610	2
Beiersdorf S.A.	Wool alcohol ointment DAB	34414610	1
Beiersdorf	Wool alcohol ointment DAB	34414610	2
Beiersdorf/Noweda	Wool alcohol ointment DAB	34424610	1
beiersdorf	Wool alcohol ointment DAB	34424610	1
Beiersdorf Ag	Wool alcohol ointment DAB	34554710	1
beiersdorf	Wool alcohol ointment DAB	40726610	1
Beiersdorf	Wool alcohol ointment DAB	50414810	1
Beiersdorf AG/ GEHE	Wool alcohol ointment DAB	41829810	1
Beiersdorf	Wool alcohol ointment DAB	50845210	1
beiersdorf	Wool alcohol ointment DAB	50815310	1
Beiersdorf	Wool alcohol ointment DAB	50955510	1
Ichthyol-Gesellschaft	Unguentum Cordes®	5052015A	1
ichthyol	Unguentum Cordes®	13A090	1
Nachprüfung	Unguentum Cordes®	12a124	1
Caelo	Unguentum Cordes®	12A211	1
Ichthyol Gesellschaft/N...	Unguentum Cordes®	13A067	1
Ichthyol -Gesellschaft	Unguentum Cordes®	13A091	1
Ichthyol / Phoenix	Unguentum Cordes®	13A091	1
Ichthyol Gesellschaft	Unguentum Cordes®	13A159	1
Ichthyol/Sanacorp	Unguentum Cordes®	13A140	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	13A180	1
Cordes, Hermani / Gehe	Unguentum Cordes®	14A112	1
Noweda/Ichthyol	Unguentum Cordes®	14A173	1
Ichthyol/Phoenix	Unguentum Cordes®	15A011	1
Ichthyol/AHP	Unguentum Cordes®	154141	1
Ichthyol	Unguentum Cordes®	14A116	1
Noweda/Ichthyol	Unguentum Cordes®	15A020	2
Alliance	Unguentum Cordes®	15A020	1
Ichthyol-Gesellschaft	Unguentum Cordes®	15A021	1
Ichthyol	Unguentum Cordes®	15A021	1
Cordes Hermann	Unguentum Cordes®	15A020	1
Ichthyol/Noweda	Unguentum Cordes®	15A023	1
Ichthyol Gesellschaft Co...	Unguentum Cordes®	15A021	1
Alliance health	Unguentum Cordes®	15A023	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/Gehe	Unguentum Cordes [®]	15A023	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	15A023	1
Phoenix	Unguentum Cordes [®]	15A027	1
Ichthol-Gesellschaft/Sa...	Unguentum Cordes [®]	15A037	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A027	1
Ichthyol/Sanova Pharma	Unguentum Cordes [®]	15A041	1
Ichthyol	Unguentum Cordes [®]	15A042	3
Caelo	Unguentum Cordes [®]	15A042	1
Ichthyol	Unguentum Cordes [®]	15A-042	1
Ichthyol/Fiebig	Unguentum Cordes [®]	15A043	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	15A043	1
Noweda / Ichthyol Gesel...	Unguentum Cordes [®]	15A060	1
Ichthyol	Unguentum Cordes [®]	15A060	2
Gehe-Ichthyol	Unguentum Cordes [®]	15A060	2
Ichthyol Gesellschaft/S...	Unguentum Cordes [®]	15A060	2
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	15A060	1
Ichthyol/Gehe	Unguentum Cordes [®]	15A060	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	15A061	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A061	1
Caelo	Unguentum Cordes [®]	15A061	2
Ichthyol/Phönix	Unguentum Cordes [®]	15A061	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A062	2
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	15A062	1
Gehe	Unguentum Cordes [®]	15A062	3
Sanacorp/Ichthyol-Gesel...	Unguentum Cordes [®]	15A062	1
Ichthyol / Sanacorp	Unguentum Cordes [®]	15A062	1
Ichthyol / Fiebig	Unguentum Cordes [®]	15A062	1
Ichthyol / Noweda	Unguentum Cordes [®]	15a073	1
Gehe-Ichthyol	Unguentum Cordes [®]	15A073	1
Caelo	Unguentum Cordes [®]	15a073	1
Ichthyol	Unguentum Cordes [®]	15A073	1
Ichthyol	Unguentum Cordes [®]	15A074	1
Ichtyol/Anzag	Unguentum Cordes [®]	15A074	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A074	1
Ichthyol/Noweda	Unguentum Cordes [®]	15A079	1
Stumpf/Ichthyol-Gesells...	Unguentum Cordes [®]	15A089	1
Ichthyol/AHC	Unguentum Cordes [®]	15A089	1
Ichthyol	Unguentum Cordes [®]	15A089	1
Ichthyol	Unguentum Cordes [®]	15A092	2
Ichthyol Gesellschaft/ ...	Unguentum Cordes [®]	15A090	1
Ichthyol über Phoenix	Unguentum Cordes [®]	15A092	1
Ichthyol	Unguentum Cordes [®]	15A093	1
Ichthyol	Unguentum Cordes [®]	15A106	1
Caelo	Unguentum Cordes [®]	15A106	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A106	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A093	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A106	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A106	1
Cordes	Unguentum Cordes [®]	15A106	1
Ichthyol	Unguentum Cordes [®]	15A107	1
phoenix Ichthyol	Unguentum Cordes [®]	15A106	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A107	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A107	1
phoenix / ichthyol	Unguentum Cordes [®]	15A108	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	15A108	1
Ichthyol/ Noweda	Unguentum Cordes [®]	15A108	1
Gehe-Ichthyol	Unguentum Cordes [®]	15A123	1
Ichthyol / Phoenix	Unguentum Cordes [®]	15A123	2
Phönix/Ichthyol	Unguentum Cordes [®]	15A123	1
Noweda/Ichthyol-Gesells...	Unguentum Cordes [®]	15A123	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/Anzag	Unguentum Cordes [®]	15A123	1
Ichthyol/Noweda	Unguentum Cordes [®]	15A123	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A124	1
Noweda/Ichthyol	Unguentum Cordes [®]	15A123	1
Ichthyol/Noweda	Unguentum Cordes [®]	15A124	1
Cordes/Phönix	Unguentum Cordes [®]	15A124	1
Caelo	Unguentum Cordes [®]	15a124	1
Noweda/Ichthyol	Unguentum Cordes [®]	15A125	3
Ichthyol Gesellschaft /...	Unguentum Cordes [®]	15A125	1
Ichthyol / Phoenix	Unguentum Cordes [®]	15A125	1
Ichthyol/Phönix	Unguentum Cordes [®]	15A125	1
Ichthyol	Unguentum Cordes [®]	15A137	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A137	1
Ichthyol-Gesellschaft/A...	Unguentum Cordes [®]	15A139	1
Ichthyol	Unguentum Cordes [®]	15A140	3
Gehe/Ichthyol	Unguentum Cordes [®]	15A141	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A141	1
Ichthyol über Phoenix	Unguentum Cordes [®]	15A142	1
Ichthyol/Fiebig	Unguentum Cordes [®]	15A140	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A142	1
Ichthyol/ Noweda	Unguentum Cordes [®]	15A148	1
ICHTHYOL-GESELLSCHAFT/P...	Unguentum Cordes [®]	15A148	1
Ichthyol- Gesellschaft	Unguentum Cordes [®]	15A149	1
Caelo	Unguentum Cordes [®]	15A149	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	15A149	2
Bombastus	Unguentum Cordes [®]	15A149	1
Phoenix	Unguentum Cordes [®]	15A149	1
Caelo	Unguentum Cordes [®]	15A150	1
Ichthyol /AHD	Unguentum Cordes [®]	15A155	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A156	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	15A159	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	15A156	1
Caelo	Unguentum Cordes [®]	15A159	1
Ichthyol / Phönix	Unguentum Cordes [®]	15A156	1
Noweda/Ichthyol	Unguentum Cordes [®]	15A191	1
Ichthyol / Gehe	Unguentum Cordes [®]	15A162	1
ichthyol-gesellschaft	Unguentum Cordes [®]	15A191	1
Ichthyol/Phoenix	Unguentum Cordes [®]	15A191	1
anzag	Unguentum Cordes [®]	15a192	1
Gehe/Ichthyol	Unguentum Cordes [®]	15A192	1
Ichtyol	Unguentum Cordes [®]	15A192	1
Noweda	Unguentum Cordes [®]	15A192	1
Ichthyol/Noweda	Unguentum Cordes [®]	15A192	2
Caelo	Unguentum Cordes [®]	15A192	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	160726UC	1
Ichthyol-Gesellschaft/...	Unguentum Cordes [®]	16A013	1
Caelo	Unguentum Cordes [®]	16A013	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A022/170117UC	1
Ichthyol über Phoenix	Unguentum Cordes [®]	16A044	1
Ichthyol/Noweda	Unguentum Cordes [®]	16A045	1
Ichthyol/ Sanacorp	Unguentum Cordes [®]	16A044	1
Ichthyol/ Krieger	Unguentum Cordes [®]	16A044	1
Ichthyol GmbH	Unguentum Cordes [®]	16A045	1
Ichthyol/N	Unguentum Cordes [®]	16A046	1
Phoenix	Unguentum Cordes [®]	16A046	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	16A047	1
Phoenix	Unguentum Cordes [®]	16A045	1
Ebert 17.08.16 KKH	Unguentum Cordes [®]	16A048	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	16A046	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A047	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	16a048	1
Caelo	Unguentum Cordes [®]	16A058	3
Ichthyol/Gehe	Unguentum Cordes [®]	16A058	1
AHD	Unguentum Cordes [®]	16A070	1
Ichthyol/Gehe	Unguentum Cordes [®]	16A083	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	16A083	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	16A058	1
Ichthyol/Anzag	Unguentum Cordes [®]	16A083	1
Gehe-Ichthyol	Unguentum Cordes [®]	16A083	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A083	2
Caelo	Unguentum Cordes [®]	16a085	1
ICHTYOL-GESELLSCHAFT Co...	Unguentum Cordes [®]	16A085	1
Ichthyol-Gesellschaft /...	Unguentum Cordes [®]	16A085	2
Cordes Hermann	Unguentum Cordes [®]	16A085	1
Gehe/Ichthyol	Unguentum Cordes [®]	16A099	3
Ichthyol/Noweda	Unguentum Cordes [®]	16A085	1
Gehe/Ichthyol	Unguentum Cordes [®]	16A085	1
Ichthyol/Noweda	Unguentum Cordes [®]	16A099	1
Ichthyol/Alliance Healt...	Unguentum Cordes [®]	16A099	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	16A099	1
Noweda/Ichthyol	Unguentum Cordes [®]	16A099	1
Ichthyol	Unguentum Cordes [®]	16A099	1
Ichthyol/Kehr	Unguentum Cordes [®]	16A101	1
Ichthyol/ Phönix	Unguentum Cordes [®]	16A101	1
Euro OTC	Unguentum Cordes [®]	16A101	1
Henry Lamotte	Unguentum Cordes [®]	16A102	1
Phoenix	Unguentum Cordes [®]	16A102	1
Sanova Pharma/Gehe	Unguentum Cordes [®]	16A102	2
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	16A102	1
Caelo	Unguentum Cordes [®]	16A119	1
Cordes/Noweda	Unguentum Cordes [®]	16A121	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A122	2
sanacorp/ichthyol	Unguentum Cordes [®]	16A121	1
ICHTHYOL-GESELLSCHAFT/P...	Unguentum Cordes [®]	16A121	1
Ichthyol / Phönix	Unguentum Cordes [®]	16A122	1
Ichthyol	Unguentum Cordes [®]	16A123	2
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	16A123	1
PHÖNIX	Unguentum Cordes [®]	16A124	1
Ichthyol	Unguentum Cordes [®]	16A124	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	16A124	1
Caelo	Unguentum Cordes [®]	16A137	2
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A137	1
Ichthyol/ Gehe	Unguentum Cordes [®]	16A139	1
Ichthyol/ Krieger	Unguentum Cordes [®]	16A139	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A139	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	16A140	1
Ichthyol	Unguentum Cordes [®]	16A140	1
Ichthyol/O.G.	Unguentum Cordes [®]	16A140	1
Ichthyol Gesellschaft/G...	Unguentum Cordes [®]	16A141	1
Cordes/Noweda	Unguentum Cordes [®]	16A140	2
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	16a141	1
sanacorp/ichthyol	Unguentum Cordes [®]	16A142	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	16A141	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	16A142	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	16A143	2
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes [®]	16A143	1
Gehe	Unguentum Cordes [®]	16A143	1
ichthyol gesellschaft	Unguentum Cordes [®]	16A143	1
Ichthyol/O.G.	Unguentum Cordes [®]	16A153	1
Ichth/Noweda	Unguentum Cordes [®]	16A144	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
ICHTHYOL-GESELLSCHAFT/P...	Unguentum Cordes®	16A144	1
Caelo	Unguentum Cordes®	16A153	2
Ichthyol/ Now	Unguentum Cordes®	16A144	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	16A155	1
Ichthyol/ Anzag	Unguentum Cordes®	16A155	1
Ichthyol-Gesellschaft	Unguentum Cordes®	16a173	1
Ichthyol- Gesellschaft	Unguentum Cordes®	16A155	1
Ichthyol/Sanacorp	Unguentum Cordes®	16A176	1
Ichthyol / Gehe	Unguentum Cordes®	16A173	1
Caelo	Unguentum Cordes®	16A176	1
Phönix / Ichthyol	Unguentum Cordes®	16A176	1
Ichtyol/Noweda	Unguentum Cordes®	16A176	1
Ichtyol	Unguentum Cordes®	16A181	1
Ichthyol/Gehe	Unguentum Cordes®	16A183	1
Bombastus	Unguentum Cordes®	16A183	1
Ichthyol	Unguentum Cordes®	16A184	1
Ichtyol	Unguentum Cordes®	16A184	1
Ichthyol/Phönix	Unguentum Cordes®	16A184	1
Sanacorp	Unguentum Cordes®	16a184	1
Cordes	Unguentum Cordes®	16A184	1
ychthyol	Unguentum Cordes®	16A184	1
Ichthyol/Sanacorp	Unguentum Cordes®	16A184	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	16A210	1
Ichthyol/Gehe	Unguentum Cordes®	16A210	1
Ichthyol/O.G.	Unguentum Cordes®	16A210	1
Ichthyol / Phönix	Unguentum Cordes®	16A210	1
Ichthyol/Noweda	Unguentum Cordes®	16A210	1
ICHTHYOL-GESELLSCHAFT/P...	Unguentum Cordes®	16A210	2
Caelo	Unguentum Cordes®	16A044	1
Ichthyol Gesellschaft /...	Unguentum Cordes®	171107UC/17A109	1
Ichthyol/Anzag	Unguentum Cordes®	1713M-02532	1
Ichthyol/OG	Unguentum Cordes®	17A003	2
Ichthyol Gesellschaft /...	Unguentum Cordes®		1
Ebert-Jacobi	Unguentum Cordes®	17A008	1
Ichthyol/Fiebig	Unguentum Cordes®	17A018	1
Ichthyol-Gesellschaft/P...	Unguentum Cordes®	17A049	1
Ichthyol/ Noweda	Unguentum Cordes®	17A002	1
Ichthyol-Gesellschaft/A...	Unguentum Cordes®	17A049	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	17A052	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	17A053	1
Ichthyol-Gesellschaft	Unguentum Cordes®	17A054	1
Ichthyol GEHE	Unguentum Cordes®	17A054	1
Ichthyol	Unguentum Cordes®	17A164	1
Ichthyol / Phönix	Unguentum Cordes®	17A164	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	17A163	1
Ichthyol/ Noweda	Unguentum Cordes®	18A002	1
ichthyol	Unguentum Cordes®		1
Fagron	Unguentum Cordes®	Ichthyol/Phoenix	1
Ichthyol-Gesellschaft	Unguentum Cordes®	31072015A	1
ICHTHYOL-GES/GEHE	Unguentum Cordes®	10A016	1
Ichthyol	Unguentum Cordes®	11A042	1
Ichthyol	Unguentum Cordes®	11A043	1
ICHTHYOL-GES/GEHE	Unguentum Cordes®	11A040	1
Ichthyol-Gesellschaft	Unguentum Cordes®	10A026	1
Ichthyol/Gehe	Unguentum Cordes®	11A043	1
Ichthyol	Unguentum Cordes®	11A120	1
Ichthyol-Gesellschaft	Unguentum Cordes®	11A085	1
Ichthyol-Gesellschaft, ...	Unguentum Cordes®	11A122	1
Caelo	Unguentum Cordes®	11A146	1
Ichthyol/Krieger	Unguentum Cordes®	11A146	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol	Unguentum Cordes [®]	11A150	1
Caelo	Unguentum Cordes [®]	11A162	1
Ichthyol-Gesellschaft, ...	Unguentum Cordes [®]	11A162	1
Ichthyol	Unguentum Cordes [®]	11A153	1
Ichthyol	Unguentum Cordes [®]	11A187	2
Sanacorp/ Sanova	Unguentum Cordes [®]	11A187	1
Noweda	Unguentum Cordes [®]	11a199	1
Ichthyol/Krieger	Unguentum Cordes [®]	11A202	1
Hersteller: Ichthyol	Unguentum Cordes [®]	11A202	1
Caelo	Unguentum Cordes [®]	11A202	1
Anzag 26.09.2012	Unguentum Cordes [®]	11A202	1
ichthyol- gesellschaft ...	Unguentum Cordes [®]	11A202	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	11A202	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	11A233	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	11A234	1
Gehe	Unguentum Cordes [®]	11A234	1
ichthyol	Unguentum Cordes [®]	12A129	1
Anzag 29.10.12	Unguentum Cordes [®]	12A007	1
Caelo	Unguentum Cordes [®]	12A021	1
Caelo	Unguentum Cordes [®]	12A056	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A041	1
Ichthyol	Unguentum Cordes [®]	12A053	2
Ichthyol/ Noweda	Unguentum Cordes [®]	12A053	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	12A053	1
Gehe	Unguentum Cordes [®]	12A054	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	12A058	1
Gehe	Unguentum Cordes [®]	12A056	2
Anzag	Unguentum Cordes [®]	12A055	1
ichthyol	Unguentum Cordes [®]	12A060	1
Caelo	Unguentum Cordes [®]	12A060	1
Ichthyol -Gesellschaft	Unguentum Cordes [®]	12A064	1
Ichthyol	Unguentum Cordes [®]	12A064	2
Gehe	Unguentum Cordes [®]	12A075	1
Ichthyol Cordes	Unguentum Cordes [®]	12A075	2
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A077	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A079	1
phönix/Ichthyol-Gesells...	Unguentum Cordes [®]	12A078	1
Caelo	Unguentum Cordes [®]	12A078	1
Gehe	Unguentum Cordes [®]	12A079	1
Phoenix	Unguentum Cordes [®]	12A087	1
Ichthyol Gesellschaft	Unguentum Cordes [®]	12A087	1
Gehe	Unguentum Cordes [®]	12A087	1
ichthyol	Unguentum Cordes [®]	12A102	1
Fagron	Unguentum Cordes [®]	12A115	1
Hersteller: Ichthyol/Gehe	Unguentum Cordes [®]	12A102	1
Stumpf/Ichthyol	Unguentum Cordes [®]	12A115	1
Caelo	Unguentum Cordes [®]	12a123	1
Ichthyol	Unguentum Cordes [®]	12A123	1
Ichthyol	Unguentum Cordes [®]	12A120	1
ichthyol,Anzag	Unguentum Cordes [®]	12A124	1
Ichthyol/Anzag	Unguentum Cordes [®]	12A124	1
Gehe	Unguentum Cordes [®]	12A124	1
Phönix	Unguentum Cordes [®]	12A129	1
Ichthyl Phönix	Unguentum Cordes [®]	12A129	1
Caelo	Unguentum Cordes [®]	12A130	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A129	1
Ichthyol	Unguentum Cordes [®]	12A139	1
Caelo	Unguentum Cordes [®]	12A131	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	12A130	1
Ichthyol Gesellschaft C...	Unguentum Cordes [®]	12A130	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol,Noweda	Unguentum Cordes®	12A139	1
Großhandel	Unguentum Cordes®	12A140	1
Ichthyol,Noweda	Unguentum Cordes®	12A142	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	12A156	1
Ichthyol	Unguentum Cordes®	12A157	1
Noweda	Unguentum Cordes®	12A184	1
Alliance	Unguentum Cordes®	12A184	1
Ichthyol-Gesellschaft C...	Unguentum Cordes®	12A184	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	12A206	1
Ichthyol/14.01.14/K	Unguentum Cordes®	13A026	1
Ichthyol Gesellschaft C...	Unguentum Cordes®	13A031	1
Gehe	Unguentum Cordes®	12A211	1
Phoenix/Ichthyol	Unguentum Cordes®	13A055	1
ichthyol	Unguentum Cordes®	13A1008	1
Ichthyol	Unguentum Cordes®	14121203	1
Ichthiol	Unguentum Cordes®	2112212532	1
Ichthyol	Unguentum Cordes®	26101201	1
Caelo	Hydrophobic base gel DAC	14102505	2
Caelo	Emulsifying hydrophobic base...	12362505	2
Caelo	Emulsifying hydrophobic base...	12303402	6
Sanacorp	Hydrophobic base gel DAC	13165201	1
Bombastus	Emulsifying hydrophobic base...	282094	1
Caelo	Hydrophilic ointment DAB	1120462	1
Alliance	Hydrophilic ointment DAB	303216009101	1
Eigen Defektur	Hydrophilic ointment DAB	316M-02539	1
Caeolo/Noweda	Hydrophilic ointment DAB	6/2019	1
Caelo	Hydrophilic ointment DAB	916I-02539	1
Caelo	Hydrophilic ointment DAB	103/AKD15	1
Caelo	Hydrophilic ointment DAB	105/AKD16	1
Caelo	Hydrophilic ointment DAB	12051401	1
Caelo	Hydrophilic ointment DAB	13117501	12
Caesar&Loretz GmbH/Anzag	Hydrophilic ointment DAB	13245301	1
Caelo	Hydrophilic ointment DAB	13245301	7
Caelo	Hydrophilic ointment DAB	13245401	10
Caelo	Hydrophilic ointment DAB	13245501	6
Caelo	Hydrophilic ointment DAB	13245701	5
Caelo	Hydrophilic ointment DAB	13245601	9
Großhandel	Hydrophilic ointment DAB	13245701	1
Caelo	Hydrophilic ointment DAB	13245901	10
Fagron	Hydrophilic ointment DAB	13H13T04006144	1
Caelo	Hydrophilic ointment DAB	14024601	4
Alliance Health 09.05.2...	Hydrophilic ointment DAB	14024701	1
Caelo	Hydrophilic ointment DAB	14024701	14
Phönix	Hydrophilic ointment DAB	14024701	1
Caelo	Hydrophilic ointment DAB	14024802	5
Caelo	Hydrophilic ointment DAB	14024803	4
AHD	Hydrophilic ointment DAB	14024803	1
Caesar & Loretz GmbH	Hydrophilic ointment DAB	14024901	2
Caelo	Hydrophilic ointment DAB	14024901	7
Caelo	Hydrophilic ointment DAB	14025001	7
Caelo	Hydrophilic ointment DAB	14123101	14
Caelo	Hydrophilic ointment DAB	14123201	4
Caelo	Hydrophilic ointment DAB	14123301	7
Caelo	Hydrophilic ointment DAB	14123401	12
Caelo	Hydrophilic ointment DAB	14123502	3
Caelo	Hydrophilic ointment DAB	14224701	4
Caelo	Hydrophilic ointment DAB	14224901	5
Caelo	Hydrophilic ointment DAB	14225001	7
Caelo	Hydrophilic ointment DAB	14225002	12
Caelo	Hydrophilic ointment DAB	1423101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic ointment DAB	14231901	8
Caelo	Hydrophilic ointment DAB	14277601	8
EuRho/ Noweda	Hydrophilic ointment DAB	14277601	1
Caelo	Hydrophilic ointment DAB	14277701	4
Caesar & Loretz GmbH	Hydrophilic ointment DAB	14277801	1
Caelo	Hydrophilic ointment DAB	14277801	7
Caelo	Hydrophilic ointment DAB	14277901	11
Caesar&Loretz GmbH/Anzag	Hydrophilic ointment DAB	14277901	1
Caesar & Loretz GmbH	Hydrophilic ointment DAB	14277901	1
Caelo	Hydrophilic ointment DAB	14278002	2
Caelo	Hydrophilic ointment DAB	14278003	3
Fagron	Hydrophilic ointment DAB	14A06-T10-009049	1
Fagron	Hydrophilic ointment DAB	14A06-T11-009054	3
Euro OTC	Hydrophilic ointment DAB	14A06-T10-009049	1
Fagron	Hydrophilic ointment DAB	14A06-T11	1
Fagron	Hydrophilic ointment DAB	14A06T11009054	1
Fagron	Hydrophilic ointment DAB	14G08-T06	1
Fagron	Hydrophilic ointment DAB	14G08-T06-013947	2
Fagron	Hydrophilic ointment DAB	14G08-T06-013948	1
Fagron	Hydrophilic ointment DAB	14G08-T07-013951	1
Fagron	Hydrophilic ointment DAB	14G08-T08-013954	4
Fagron	Hydrophilic ointment DAB	14G08-To8-013954	1
Caelo	Hydrophilic ointment DAB	15016502	7
Caelo	Hydrophilic ointment DAB	15016501	1
Caelo	Hydrophilic ointment DAB	15016503	4
Caelo	Hydrophilic ointment DAB	15016601	11
Caelo	Hydrophilic ointment DAB	15016701	5
Caelo	Hydrophilic ointment DAB	15016602	2
Caelo	Hydrophilic ointment DAB	15016801	8
Caelo	Hydrophilic ointment DAB	15016901	10
Caelo	Hydrophilic ointment DAB	15042014B	1
Caelo	Hydrophilic ointment DAB	15097714	1
Caelo	Hydrophilic ointment DAB	15113602	1
Caelo	Hydrophilic ointment DAB	15221501	5
Caelo	Hydrophilic ointment DAB	15221601	10
Noweda 06.04.2016	Hydrophilic ointment DAB	15221601	1
Cealo/Noweda	Hydrophilic ointment DAB	15221601	1
Caelo	Hydrophilic ointment DAB	15221602	6
Caelo	Hydrophilic ointment DAB	15221701	19
Caesar & Loretz GmbH	Hydrophilic ointment DAB	15221701	2
Caelo	Hydrophilic ointment DAB	15221801	2
Caelo	Hydrophilic ointment DAB	15221901	2
Caelo	Hydrophilic ointment DAB	15221902	2
Caelo	Hydrophilic ointment DAB	15400701	1
Caelo	Hydrophilic ointment DAB	15416401	12
Caelo	Hydrophilic ointment DAB	15416502	10
Fagron	Hydrophilic ointment DAB	15B19-T01	1
Fagron	Hydrophilic ointment DAB	15B19-T01-019411	2
Fagron	Hydrophilic ointment DAB	15B19-T01-019409	1
sanacorp	Hydrophilic ointment DAB	15F09-T02023370	1
Fagron	Hydrophilic ointment DAB	15F09-T02-023370	1
Fagron	Hydrophilic ointment DAB	15F24-T05-023887	1
Fagron	Hydrophilic ointment DAB	15K24-T01-027590	1
Caelo	Hydrophilic ointment DAB	16008701	6
Caesar & Loretz GmbH Hi...	Hydrophilic ointment DAB	16008701	2
phönix	Hydrophilic ointment DAB	16008701	1
Caelo	Hydrophilic ointment DAB	16008702	7
Caelo	Hydrophilic ointment DAB	16008901	8
Caelo	Hydrophilic ointment DAB	16009001	7
Ebert & Jacobi 03.08.16	Hydrophilic ointment DAB	160091001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic ointment DAB	16009101	17
Caesar&Loretz GmbH/Anzag	Hydrophilic ointment DAB	16009101	1
Celo	Hydrophilic ointment DAB	16009101	1
Caelo	Hydrophilic ointment DAB	16009201	8
Anzag/Ichtyol-Gesellsch. ...	Hydrophilic ointment DAB	16009201	1
Caelo	Hydrophilic ointment DAB	16009301	5
Caelo	Hydrophilic ointment DAB	16009401	2
Caelo	Hydrophilic ointment DAB	160408UE	1
Caelo	Hydrophilic ointment DAB	1614U2539	1
Caelo	Hydrophilic ointment DAB	16166701	10
Caesar & Loretz	Hydrophilic ointment DAB	16166701	1
phönix	Hydrophilic ointment DAB	16166801	1
Caelo	Hydrophilic ointment DAB	16166801	22
PhönixWE:2.02.2017 EK...	Hydrophilic ointment DAB	16166801	1
Bombastus	Hydrophilic ointment DAB	16166801	1
Caelo	Hydrophilic ointment DAB	16166901	4
Phönix	Hydrophilic ointment DAB	16166901	1
Caelo	Hydrophilic ointment DAB	16166902	6
Caelo	Hydrophilic ointment DAB	16167001	6
Sanacorp 18.03.17 Ek:11...	Hydrophilic ointment DAB	16167001	1
Gehe	Hydrophilic ointment DAB	16167001	1
Caeolo/Noweda	Hydrophilic ointment DAB	16167102	1
Caelo	Hydrophilic ointment DAB	1617102	2
Caelo	Hydrophilic ointment DAB	16282401	2
phönix	Hydrophilic ointment DAB	16282401	1
Caesar und Loretz GmbH	Hydrophilic ointment DAB	16282401	2
Caelo	Hydrophilic ointment DAB	16282402	5
Caelo	Hydrophilic ointment DAB	16282501	7
Phönix	Hydrophilic ointment DAB	16282702	1
Caelo	Hydrophilic ointment DAB	16282801	2
Caelo	Hydrophilic ointment DAB	1628401	1
Fagron	Hydrophilic ointment DAB	16G18-T04-035261	3
Caelo	Hydrophilic ointment DAB	16K09-T03-03854	1
Caelo	Hydrophilic ointment DAB	17040201	6
Caleo/Noweda	Hydrophilic ointment DAB	17040201	1
Caelo	Hydrophilic ointment DAB	17040202	2
phönix	Hydrophilic ointment DAB	17040301	1
Caelo	Hydrophilic ointment DAB	17040301	2
Wepa/ Phoenix	Hydrophilic ointment DAB	17040301	1
Allianz Healthcare	Hydrophilic ointment DAB	17040301	1
Caelo	Hydrophilic ointment DAB	17040501	2
Caelo	Hydrophilic ointment DAB	17040601	5
Caelo	Hydrophilic ointment DAB	17083002	2
Caelo	Hydrophilic ointment DAB	17083201	4
Caelo	Hydrophilic ointment DAB	17083004	1
Caelo	Hydrophilic ointment DAB	17083301	1
Caelo	Hydrophilic ointment DAB	170905UE/16167001	1
Eigen Defektur	Hydrophilic ointment DAB	1716E-02539	1
Caelo	Hydrophilic ointment DAB	17185501	3
Aponorm/Wepa	Hydrophilic ointment DAB	17185501	1
Caelo	Hydrophilic ointment DAB	17290603	1
Caelo	Hydrophilic ointment DAB	17290701	3
Noweda	Hydrophilic ointment DAB	17290702	1
Caelo	Hydrophilic ointment DAB	17290703	1
Caelo	Hydrophilic ointment DAB	17290702	1
Caelo	Hydrophilic ointment DAB	17290801	1
Fagron	Hydrophilic ointment DAB	17I11-T03-048188	1
Henry Lamotte/Anzag	Hydrophilic ointment DAB	2215M-02539	1
Eigen/Defektur	Hydrophilic ointment DAB	2516A-02539	1
Bombastus	Hydrophilic ointment DAB	279002	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hydrophilic ointment DAB	281346	2
Bombastus	Hydrophilic ointment DAB	281248	1
Bombastus	Hydrophilic ointment DAB	281347	3
Bombastus	Hydrophilic ointment DAB	281348	3
Caelo	Hydrophilic ointment DAB	281349	2
Bombastus	Hydrophilic ointment DAB	281351	6
Bombastus	Hydrophilic ointment DAB	284746	2
Bombastus	Hydrophilic ointment DAB	284747	1
Bombastus	Hydrophilic ointment DAB	286544	1
Caelo	Hydrophilic ointment DAB	286544	5
Bombastus	Hydrophilic ointment DAB	286545	9
Caelo	Hydrophilic ointment DAB	286545	1
Bombastus	Hydrophilic ointment DAB	288366	1
Bombastus	Hydrophilic ointment DAB	288367	4
Bombastus	Hydrophilic ointment DAB	288368	1
Bombastus	Hydrophilic ointment DAB	288369	6
Euro OTC	Hydrophilic ointment DAB	288369	1
Bombastus	Hydrophilic ointment DAB	290597	2
Bombastus	Hydrophilic ointment DAB	290599	4
Bombastus	Hydrophilic ointment DAB	290600	3
Bombastus	Hydrophilic ointment DAB	290601	9
Bombastus	Hydrophilic ointment DAB	293261	4
Bombastus	Hydrophilic ointment DAB	293262	4
Bombastus	Hydrophilic ointment DAB	293263	2
Anzag	Hydrophilic ointment DAB	293264	1
Bombastus	Hydrophilic ointment DAB	293264	6
Bombastus	Hydrophilic ointment DAB	294448	4
Caelo	Hydrophilic ointment DAB	294449	2
Bombastus	Hydrophilic ointment DAB	294451	3
Bombastus	Hydrophilic ointment DAB	296866	3
Bombastus	Hydrophilic ointment DAB	296865	3
Bombastus	Hydrophilic ointment DAB	296868	6
Bombastus	Hydrophilic ointment DAB	296867	1
Defektur / Wolf-Apotheke	Hydrophilic ointment DAB	299390	3
Bombastus	Hydrophilic ointment DAB	299390	2
Bombastus	Hydrophilic ointment DAB	299392	3
Bombastus	Hydrophilic ointment DAB	299391	2
Caelo	Hydrophilic ointment DAB	299392	1
Bombastus	Hydrophilic ointment DAB	301096	3
Bombastus	Hydrophilic ointment DAB	301094	1
Bombastus	Hydrophilic ointment DAB	301095	1
Bombastus	Hydrophilic ointment DAB	301097	2
Caelo	Hydrophilic ointment DAB	301097	1
Bombastus	Hydrophilic ointment DAB	303000	1
Bombastus	Hydrophilic ointment DAB	302999	1
Bombastus	Hydrophilic ointment DAB	302998	1
Bombastus	Hydrophilic ointment DAB	305394	2
Bombastus	Hydrophilic ointment DAB	305395	4
Caelo	Hydrophilic ointment DAB	305396	1
Bombastus	Hydrophilic ointment DAB	307092	1
Bombastus	Hydrophilic ointment DAB	307094	2
Bombastus	Hydrophilic ointment DAB	311063	1
Defekur	Hydrophilic ointment DAB	3216M-02539	1
Eigen/Defektur	Hydrophilic ointment DAB	3415I-02539	1
Eigen/Defektur	Hydrophilic ointment DAB	4215E-02539	1
Eigen/Defektur	Hydrophilic ointment DAB	5015E-02539	1
Anzag 17.03.2014 EK 13...	Hydrophilic ointment DAB	5448000	1
lamotte	Hydrophilic ointment DAB	5448000	1
Lamotte/Alliance	Hydrophilic ointment DAB	5651600	1
lamotte	Hydrophilic ointment DAB	5517900	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fiebig/ Lamotte	Hydrophilic ointment DAB	5860400	1
Lamotte//AHD	Hydrophilic ointment DAB	5860500	1
Lamotte/ANZAG	Hydrophilic ointment DAB	5860400	1
Henry Lamotte Oils Gmbh	Hydrophilic ointment DAB	5862400	1
Henry Lamotte/Noweda	Hydrophilic ointment DAB	6001368014	2
Sanacorp	Hydrophilic ointment DAB	6/AKD7	1
Caelo	Hydrophilic ointment DAB		4
Sanacorp	Hydrophilic ointment DAB		1
Caelo	Hydrophilic ointment DAB	11133502	1
Caelo	Hydrophilic ointment DAB	11211601	2
Caelo	Hydrophilic ointment DAB	11212101	5
Caelo	Hydrophilic ointment DAB	11133701	1
Caelo	Hydrophilic ointment DAB	11212201	2
Gehe	Hydrophilic ointment DAB	11L13-T01	2
Caelo	Hydrophilic ointment DAB	12014601	1
Gehe	Hydrophilic ointment DAB	12014601	1
Caelo	Hydrophilic ointment DAB	12014602	3
Caelo	Hydrophilic ointment DAB	2803212539	1
Caelo	Hydrophilic ointment DAB	12014801	1
Caelo	Hydrophilic ointment DAB	12014701	3
Bombastus	Hydrophilic ointment DAB	269532	1
Bombastus	Hydrophilic ointment DAB	269536	1
Bombastus	Hydrophilic ointment DAB	269537	1
Bombastus	Hydrophilic ointment DAB	269538	3
Stumpf	Hydrophilic ointment DAB	269539	1
Bombastus	Hydrophilic ointment DAB	269539	4
Bombastus	Hydrophilic ointment DAB	274004	1
Caelo	Hydrophilic ointment DAB	103312539	1
Stumpf	Hydrophilic ointment DAB	274004	1
Caelo	Hydrophilic ointment DAB	10025601	1
Caelo	Hydrophilic ointment DAB	11015402	1
Caelo	Hydrophilic ointment DAB	11133601	3
Caelo	Hydrophilic ointment DAB	70913UngEm	1
Caelo	Hydrophilic ointment DAB	10082901	1
Euro OTC	Hydrophilic ointment DAB	11211601	2
Caelo	Hydrophilic ointment DAB	11212001	1
Caelo	Hydrophilic ointment DAB	11335101	2
Noweda	Hydrophilic ointment DAB	12014601	1
Caelo	Hydrophilic ointment DAB	12055801	3
Caelo	Hydrophilic ointment DAB	12015001	3
Caelo	Hydrophilic ointment DAB	12055901	10
Caelo	Hydrophilic ointment DAB	12056001	1
Caelo	Hydrophilic ointment DAB	12056002	4
Caesar & Loretz GmbH, 4...	Hydrophilic ointment DAB	12056002	1
Fiebig 08.10.12 ek: 10,...	Hydrophilic ointment DAB	12055901	1
Caelo	Hydrophilic ointment DAB	12086702	3
Caelo	Hydrophilic ointment DAB	12204802	8
Noweda	Hydrophilic ointment DAB	12204802	2
Caelo	Hydrophilic ointment DAB	12204803	1
Caelo	Hydrophilic ointment DAB	12205001	6
Caelo	Hydrophilic ointment DAB	12204901	2
Sanacorp	Hydrophilic ointment DAB	12205001	1
Caelo	Hydrophilic ointment DAB	12205101	4
Großhandel	Hydrophilic ointment DAB	12205101	1
Caelo	Hydrophilic ointment DAB	12205201	7
Gehe	Hydrophilic ointment DAB	12205201	1
Noweda	Hydrophilic ointment DAB	12205201	1
Caelo	Hydrophilic ointment DAB	12311102	5
Alliance Healthcare	Hydrophilic ointment DAB	12311102	1
Phönix	Hydrophilic ointment DAB	12311301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic ointment DAB	12311301	7
Caelo	Hydrophilic ointment DAB	12311401	12
Krieger	Hydrophilic ointment DAB	12311401	1
Caelo	Hydrophilic ointment DAB	12311501	16
Doroderm GmbH/ Noweda	Hydrophilic ointment DAB	12311501	1
Caelo	Hydrophilic ointment DAB	12311502	2
Caesar&lorenz,Phoenix	Hydrophilic ointment DAB	12311501	1
Noweda	Hydrophilic ointment DAB	12311502	1
Caelo	Hydrophilic ointment DAB	12311701	7
Caelo	Hydrophilic ointment DAB	12311801	3
Noweda	Hydrophilic ointment DAB	12311801	1
Fagron	Hydrophilic ointment DAB	12A12-T03	1
Fagron	Hydrophilic ointment DAB	12C09-T01	2
Phönix Gotha	Hydrophilic ointment DAB	12G18-T10	1
Fagron	Hydrophilic ointment DAB	12J31-T04	4
Caelo	Hydrophilic ointment DAB	12J31T04	1
Caelo	Hydrophilic ointment DAB	13013901	2
Caelo	Hydrophilic ointment DAB	13014001	4
Caelo	Hydrophilic ointment DAB	13014002	9
Gehe 07.05.2013	Hydrophilic ointment DAB	13014002	1
Caesar & Loretz GmbH	Hydrophilic ointment DAB	13014002	1
Caelo	Hydrophilic ointment DAB	13014003	5
Caelo	Hydrophilic ointment DAB	13014101	6
gehe 11.06.2013	Hydrophilic ointment DAB	13014101	1
Euro OTC	Hydrophilic ointment DAB	13117102	1
Caelo	Hydrophilic ointment DAB	13117102	10
Sanacorp	Hydrophilic ointment DAB	13117102	3
Caelo	Hydrophilic ointment DAB	13117103	4
Caelo	Hydrophilic ointment DAB	13117201	2
Caelo	Hydrophilic ointment DAB	13117301	20
Gehe	Hydrophilic ointment DAB	13117301	2
Anzag 02.10.13	Hydrophilic ointment DAB	13117301	1
Caelo	Hydrophilic ointment DAB	13117401	7
Caesar & Loretz GmbH, 4...	Hydrophilic ointment DAB	13117501	1
Großhandel	Hydrophilic ointment DAB	13117501	1
Alliance Healthcare	Hydrophilic ointment DAB	13117501	1
Caelo	Hydrophilic ointment DAB	13117601	6
Caelo	Hydrophilic ointment DAB	13245302	7
Caelo	Hydrophilic ointment DAB	13245303	1
Caesar&lorenz,Phoenix	Hydrophilic ointment DAB	13245401	1
Caesar&Loretz GmbH/Anzag	Hydrophilic ointment DAB	13245401	1
Caesar & Loretz GmbH	Hydrophilic ointment DAB	13245401	1
Caelo	Hydrophilic ointment DAB	13245801	4
Caelo	Hydrophilic ointment DAB	13246002	4
Fagron	Hydrophilic ointment DAB	13B04-T18-001057	3
Caelo	Hydrophilic ointment DAB	13B04-T18-001061	2
Fagron	Hydrophilic ointment DAB	13E10-T06-003561	1
Fagron	Hydrophilic ointment DAB	13E10-T06-003563	1
Fagron	Hydrophilic ointment DAB	13H13-T04-006144	1
Bombastus	Hydrophilic ointment DAB	14025001	2
Caelo	Hydrophilic ointment DAB	15056002	1
Caelo	Hydrophilic ointment DAB	2109212539	1
Caelo	Hydrophilic ointment DAB	25051307	1
Caelo	Hydrophilic ointment DAB	20221110-1	1
Caelo	Hydrophilic ointment DAB	24011315	1
Caelo	Hydrophilic ointment DAB	27121301	1
Bombastus	Hydrophilic ointment DAB	276671	5
Gehe	Hydrophilic ointment DAB	276671	2
Bombastus	Hydrophilic ointment DAB	279001	2
Bombastus	Hydrophilic ointment DAB	281349	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophilic ointment DAB	281351	1
Bombastus	Hydrophilic ointment DAB	281350	1
Caelo	Hydrophilic ointment DAB	3613Q2539	1
Caelo	Hydrophilic ointment DAB	4912E-02539	1
Caelo	Hydrophilic ointment DAB	3962/V04/16	1
Lamotte	Hydrophilic ointment DAB	518400	1
Caelo	Hydrophilic ointment DAB	5099200	1
Phönix	Hydrophilic ointment DAB	5321400	1
DAB 2008 Henry Lamotte ...	Hydrophilic ointment DAB	5243000	1
Henry Lamotte/Anzag	Hydrophilic ointment DAB	5321400	1
Lamotte	Hydrophilic ointment DAB	5321400	1
Caelo	Hydrophilic ointment DAB	5448000	1
Bombastus	Hydrophilic ointment DAB	6861252	1
Lamotte/Alliance	Hydrophilic ointment DAB	5517900	1
Caelo	Hydrophilic ointment DAB	80127028	1
Lamotte / AHC	Hydrophilic ointment DAB	5448000	1
Sanacorp	Hydrophilic ointment DAB	AKD32	1
Phoenix*lamotte	Hydrophilic ointment DAB	5517900	1
Caelo	Hydrophilic ointment DAB	Test	1
Fagron	Hydrophilic ointment DAB		1
Caelo	Emulsifying hydrophobic base...		1
Gehe Kassel	Hydrophobic base gel DAC	271904	1
Caelo	Hydrophilic ointment DAB	U-020413-1	1
Bombastus	Hydrophobic base gel DAC	10071402	1
Fagron	Hydrophobic base gel DAC	10L14-T54-006109	4
Caelo	Hydrophobic base gel DAC	64/2014	1
Fagron	Emulsifying hydrophobic base...	10L14-T54-003418	1
Fagron	Hydrophobic base gel DAC	10L14-T54-009128	2
Caelo	Hydrophobic base gel DAC	10L14-T54-009128	1
Caelo	Hydrophobic base gel DAC	13341701	1
Gehe 23.05.2014	Hydrophobic base gel DAC	14044601	1
Phoenix	Emulsifying hydrophobic base...	14087502	1
Caelo	Hydrophobic base gel DAC	14304603	3
Caelo	Hydrophobic base gel DAC	14304624	1
Caelo	Hydrophobic base gel DAC	14304604	1
Caelo	Hydrophobic base gel DAC	15055702	2
Noweda	Hydrophobic base gel DAC	15055705	1
Caelo	Emulsifying hydrophobic base...	15064001	3
Caelo	Emulsifying hydrophobic base...	15064201	3
Geilenkirchen	Hydrophobic base gel DAC	15064201	1
Caelo	Hydrophobic base gel DAC	15149602	1
Noweda	Hydrophobic base gel DAC	15149604	1
Caelo	Emulsifying hydrophobic base...	15205101	1
Caelo	Emulsifying hydrophobic base...	15205102	1
Noweda	Hydrophobic base gel DAC	15240805	2
Caelo	Hydrophobic base gel DAC	15240802	2
Fagron	Hydrophobic base gel DAC	15240805	1
Caelo	Hydrophobic base gel DAC	15240807	3
Caelo	Emulsifying hydrophobic base...	15269102	1
Noweda	Hydrophobic base gel DAC	15331102	1
Caelo	Hydrophobic base gel DAC	15383005	1
Fagron	Hydrophobic base gel DAC	15436101	3
Caelo	Hydrophobic base gel DAC	15436102	4
Euro OTC	Hydrophobic base gel DAC	15436102	1
Caelo	Hydrophobic base gel DAC	15436105	2
Caesar/Noweda	Hydrophobic base gel DAC	15436103	1
Caelo	Emulsifying hydrophobic base...	15964101	1
Fagron	Hydrophobic base gel DAC	15J02-T09-023890	1
Fagron	Hydrophobic base gel DAC	15L10-T04-028696	2
Bombastus	Hydrophobic base gel DAC	160607HB	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrophobic base gel DAC	16151101	2
Caelo	Hydrophobic base gel DAC	16151105	5
Caelo	Emulsifying hydrophobic base...	16151102	1
Gehe	Hydrophobic base gel DAC	16151105	1
Celo	Hydrophobic base gel DAC	16172301	1
Caelo	Hydrophobic base gel DAC	16172301	1
Caelo	Emulsifying hydrophobic base...	16226002	3
Caelo	Hydrophobic base gel DAC	16240105	3
Gehe	Hydrophobic base gel DAC	16240104	1
Caelo	Hydrophobic base gel DAC	16240110	1
Caelo	Hydrophobic base gel DAC	16240111	3
Caelo	Hydrophobic base gel DAC	16240114	2
Fagron	Hydrophobic base gel DAC		1
Caelo	Emulsifying hydrophobic base...	17032102	1
Caelo	Hydrophobic base gel DAC	17032201	1
Caelo	Hydrophobic base gel DAC	17043605	2
Gehe	Hydrophobic base gel DAC	17113702	1
Caelo	Emulsifying hydrophobic base...	17113702	1
Caelo	Hydrophobic base gel DAC	17196704	1
Caelo	Hydrophobic base gel DAC	17196713	1
Caelo	Hydrophobic base gel DAC	17196703	1
Bombastus	Hydrophobic base gel DAC	24031504	1
Phoenix 02.05.2014	Emulsifying hydrophobic base...	282094	1
Bombastus	Hydrophobic base gel DAC	283998	1
Bombastus	Hydrophobic base gel DAC	286013	1
GEHE Kassel	Hydrophobic base gel DAC	289214	1
Bombastus	Hydrophobic base gel DAC	292717	1
Bombastus	Hydrophobic base gel DAC	297433	1
Bombastus	Hydrophobic base gel DAC	299263	1
Bombastus	Hydrophobic base gel DAC	300895	1
Bombastus	Hydrophobic base gel DAC	302145	1
Caelo	Hydrophobic base gel DAC	302145	1
Bombastus	Hydrophobic base gel DAC	302145	2
Bombastus	Hydrophobic base gel DAC	304035	1
Sanacorp	Hydrophobic base gel DAC	308440	1
Bombastus	Hydrophobic base gel DAC	271904	2
Bombastus	Emulsifying hydrophobic base...	271106	2
Bombastus	Emulsifying hydrophobic base...	262696	1
Bombastus	Emulsifying hydrophobic base...	255106	1
Anzag 05.11.12	Emulsifying hydrophobic base...	271106	1
Fagron	Emulsifying hydrophobic base...	10L14-T54-111713	1
Caelo	Emulsifying hydrophobic base...	11016804	1
Caelo	Emulsifying hydrophobic base...	11105002	1
Bombastus	Emulsifying hydrophobic base...	11041301	1
Bombastus	Emulsifying hydrophobic base...	11071308	1
Fagron	Emulsifying hydrophobic base...	1108-T01	1
Caelo	Emulsifying hydrophobic base...	11209302	1
Fagron	Emulsifying hydrophobic base...	11I08-T01	2
Caelo	Emulsifying hydrophobic base...	11301605	1
Fagron	Emulsifying hydrophobic base...	11G25-T01	1
Caelo	Emulsifying hydrophobic base...	12078001	2
Caelo	Emulsifying hydrophobic base...	12160201	4
Noweda	Emulsifying hydrophobic base...	12160201	1
Caelo	Emulsifying hydrophobic base...	12182303	3
phoenix	Emulsifying hydrophobic base...	12182301	1
Caelo	Emulsifying hydrophobic base...	12229001	4
anzag,04.12.12	Emulsifying hydrophobic base...	12229006	1
phönix, 21.01.2013	Emulsifying hydrophobic base...	12303402	1
Caelo	Emulsifying hydrophobic base...	12362504	1
Fagron	Emulsifying hydrophobic base...	12F29-T01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Hydrophobic base gel DAC	12F29-T01	1
Gehe	Emulsifying hydrophobic base...	13066802	3
Caelo	Hydrophobic base gel DAC	13165205	2
Gehe	Emulsifying hydrophobic base...	13177601	1
Sanacorp	Emulsifying hydrophobic base...	13177601	1
Caelo	Emulsifying hydrophobic base...	13177602	1
Caelo	Hydrophobic base gel DAC	13295701	1
Caelo	Emulsifying hydrophobic base...	13341701	2
Caelo	Emulsifying hydrophobic base...	13341703	3
Noweda	Emulsifying hydrophobic base...	13361004	1
Bombastus	Emulsifying hydrophobic base...	22111307	1
Bombastus	Emulsifying hydrophobic base...	277443	1
Caelo	Emulsifying hydrophobic base...	28051301	1
Apomix	Eye ointment, emulsifying DAC	3A180213	1
Caelo	Paraffin, white	13268422	5
Euro OTC	Wool alcohol ointment DAB	1211017-01	3
Caelo	Wool alcohol ointment DAB	12340101	6
Caelo	Wool alcohol ointment DAB	15035502	13
Caelo	Wool alcohol ointment DAB	15225401	2
Caelo	Wool alcohol ointment DAB	15228201	8
Beiersdorf AG	Wool alcohol ointment DAB	15362002	1
Caelo	Wool alcohol ointment DAB	17104902	3
Caelo	Wool alcohol ointment DAB	16198502	12
Bombastus	Wool alcohol ointment DAB	296437	2
Bombastus	Wool alcohol ointment DAB	295218	3
Bombastus	Wool alcohol ointment DAB	301202	4
Beiersdorf/Phönix	Wool alcohol ointment DAB	60429110	1
Anzag/Beiersdorf	Wool alcohol ointment DAB	50815310	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	63052310	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	61149910	1
Bombastus	Wool alcohol ointment DAB	265368	1
Bombastus	Wool alcohol ointment DAB	270643	4
Caelo	Wool alcohol ointment DAB	11335002	6
Holdermann/Beiersdorf	Wool alcohol ointment DAB	22532710	1
Caelo	Wool alcohol ointment DAB	24243910	1
Beiersdorf/N/22.03.2013...	Wool alcohol ointment DAB	24635410	1
Apomix/Sanacorp 13.05.2...	Eye ointment, emulsifying DAC	3270115	1
AHD /apomix	Eye ointment, emulsifying DAC	3B060217	1
Apomix/ Phönix	Eye ointment, emulsifying DAC	3C200114	2
apomix/Phönix	Eye ointment, emulsifying DAC	4C060217	1
Apomix/ Phönix	Eye ointment, emulsifying DAC	7A010616	1
Apomix / Phönix	Eye ointment, emulsifying DAC	7B030314	1
Apomix /Noweda	Eye ointment, emulsifying DAC	9B050916	1
Apomix /Noweda	Eye ointment, emulsifying DAC	10C041016	1
Apomix/Kehr	Eye ointment, emulsifying DAC	11B011116	1
apomix	Eye ointment, emulsifying DAC	12C010814	1
Caelo	Eye ointment, emulsifying DAC	14C030815	1
PKH/ Noweda	Eye ointment, emulsifying DAC	12E010814	1
Apomix / Phönix	Eye ointment, emulsifying DAC	15A170914	1
Apomix/Kehr	Eye ointment, emulsifying DAC	15B030815	1
Apomix/ Phönix	Eye ointment, emulsifying DAC	13A010715	1
Apomix/Noweda	Eye ointment, emulsifying DAC	15C030815	1
apomix	Eye ointment, emulsifying DAC	15D030815	1
Euro OTC	Eye ointment, emulsifying DAC	18A021115	1
Apomix/ Phönix	Eye ointment, emulsifying DAC	19A021214	1
AMH Niemann/E	Eye ointment, emulsifying DAC	1C080114	1
Apomix/Kehr	Eye ointment, emulsifying DAC	19B021115	1
PKH Halle/ Phoenix	Eye ointment, emulsifying DAC	3A180213	1
Phönix Gotha	Eye ointment, emulsifying DAC	5B220413	1
PKH Gmbh/Phönix	Eye ointment, emulsifying DAC	7B210513	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PKH Halle/ Phoenix	Eye ointment, emulsifying DAC	11D170713	1
apomix/Phönix	Eye ointment, emulsifying DAC	11A170713	1
PKH GmbH/ Phönix	Eye ointment, emulsifying DAC	12C300813	1
Phönix Gotha	Eye ointment, emulsifying DAC	16B061212	1
Apomix /	Eye ointment, emulsifying DAC	17B111113	1
Apomix / Phönix	Eye ointment, emulsifying DAC	17C061212	2
Caelo	Paraffin, liquid		2
Bombastus	Paraffin, liquid	260652	1
Bombastus	Paraffin, liquid	263350	1
Bombastus	Paraffin, liquid	266032	1
Phoenix, 04.08.2012	Paraffin, liquid	269213	2
Bombastus	Paraffin, liquid	272590	1
Fagron	Paraffin, liquid	941435	1
Henry Lamotte	Paraffin, liquid	122587	1
Noweda	Paraffin, liquid	122875	1
Gehe	Paraffin, liquid	122875	1
Henry Lamotte GmbH, Mer...	Paraffin, liquid	122875	1
Henry Lamotte Oils GmbH	Paraffin, liquid	122875	1
Henry Lamotte	Paraffin, liquid	122875	2
Lamotte	Paraffin, liquid	123075	1
H.Lamotte/Noweda	Paraffin, liquid	123075	1
Noweda / Henry Lamotte	Paraffin, liquid	123075	2
Henry Lamotte	Paraffin, liquid	123075	1
Phönix	Paraffin, liquid	123075	1
Henry Lamotte Oils GmbH	Paraffin, liquid	123075	1
Henry Lamotte oils GmbH	Paraffin, liquid	123075	2
Henry Lamotte	Paraffin, liquid	123393	1
Phönix	Paraffin, liquid	123393	2
Henry Lamotte Oels GmbH	Paraffin, liquid	123393	1
Henry Lamotte	Paraffin, liquid	123743	11
Henry Lamotte oils GmbH	Paraffin, liquid	123743	8
Lamotte	Paraffin, liquid	123743	1
Henry Lamotte Oils GmbH...	Paraffin, liquid	123743	3
Henry Lamotte / Noweda	Paraffin, liquid	123743	4
Henry Lamotte GmbH, Mer...	Paraffin, liquid	123743	1
Noweda	Paraffin, liquid	123743	1
Henry Lamotte Oils / No...	Paraffin, liquid	123743	2
Oils/ Phönix	Paraffin, liquid	123743	2
Henry Lamotte Oils GmbH...	Paraffin, liquid	123743	1
Lamotte/Anzag	Paraffin, liquid	123743	1
Caelo	Paraffin, liquid	123743	1
Lamotte / Alliance Heal...	Paraffin, liquid	125544	1
Caelo	Paraffin, liquid	2111202	1
Caelo	Paraffin, liquid	214U-02050	1
Caelo	Paraffin, liquid	4091515	1
Eigenherstellung	Paraffin, liquid	60614HJK01	1
Bombastus	Paraffin, liquid	6861849	2
Caelo	Paraffin, liquid	715I-02050	1
Caelo	Paraffin, liquid	9320905	1
henry lamotte oils/Sana...	Paraffin, liquid	1000002995	1
Fagron	Paraffin, liquid	1004C014	1
Eigenherstellung	Paraffin, liquid	100615HJK01	1
Caelo	Paraffin, liquid	10062014A	1
Caelo	Paraffin, liquid	10071202	1
Caelo	Paraffin, liquid	10071403	1
Caelo	Paraffin, liquid	10162002	1
Caelo	Paraffin, liquid	10162009	1
Caelo	Paraffin, liquid	10364703	1
Caelo	Paraffin, liquid	11009510	1
Euro OTC	Paraffin, liquid	1106009-01	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Paraffin, liquid	1108033-01	1
Caelo	Paraffin, liquid	11071302	1
Caelo	Paraffin, liquid	11092703	1
Fagron	Paraffin, liquid	110816-2	1
Euro OTC	Paraffin, liquid	1110014-01	3
Hersteller: EuRho	Paraffin, liquid	1110014-01	1
Europharm/Otto	Paraffin, liquid	1110014-01	1
Caelo	Paraffin, liquid	1114U02050	1
Caelo	Paraffin, liquid	11207702	3
Eigenherstellung	Paraffin, liquid	11214HAS01	1
Caelo	Paraffin, liquid	1126261	1
Caelo	Paraffin, liquid	11263606	4
Caelo	Paraffin, liquid	11354102	7
Gehe	Paraffin, liquid	11354102	1
Noweda	Paraffin, liquid	11354102	1
Caelo	Paraffin, liquid	1135412	1
Caelo	Paraffin, liquid	113544102	1
Fagron	Paraffin, liquid	11co8-no7	1
Caelo	Paraffin, liquid	11H04-N02	1
Kehr	Paraffin, liquid	11H04-N02	1
Fagron	Paraffin, liquid	11L19-N01	1
Fagron	Paraffin, liquid	11f30-n03	1
Euro OTC	Paraffin, liquid	1203042-01	4
EuRho	Paraffin, liquid	1203042-01	1
Phoenix	Paraffin, liquid	12042202	1
Caelo	Paraffin, liquid	12042202	19
Sanacorp	Paraffin, liquid	12042202	3
Anzag 02.01.2013	Paraffin, liquid	12042202	1
Caelo	Paraffin, liquid	12042205	4
Sanacorp	Paraffin, liquid	12042205	2
Caelo	Paraffin, liquid	1205033-01	1
Gehe	Paraffin, liquid	1207022-01	2
EuRho	Paraffin, liquid	1207022-01	1
Euro OTC	Paraffin, liquid	1207022-01	2
Euro OTC	Paraffin, liquid	1208021-01	3
Noweda	Paraffin, liquid	1208021-01	3
Gehe	Paraffin, liquid	1208021-01	1
Euro OTC	Paraffin, liquid	120802101	1
Caelo	Paraffin, liquid	12102101	11
Caelo	Paraffin, liquid	12109603	2
Euro OTC	Paraffin, liquid	1211026-01	7
Euro OTC	Paraffin, liquid	1211027-01	9
Gehe	Paraffin, liquid	1211026-01	1
Phönix/Euro	Paraffin, liquid	1211027-01	1
Europharm/Otto	Paraffin, liquid	1211027-01	1
Noweda	Paraffin, liquid	1211027-01	1
Euro OTC	Paraffin, liquid	121102701	1
Caelo	Paraffin, liquid	12227601	2
Caelo	Paraffin, liquid	12227602	1
Caelo	Paraffin, liquid	12227609	6
Gehe 16.04.2013	Paraffin, liquid	12227609	1
Caelo	Paraffin, liquid	12227610	1
Caelo	Paraffin, liquid	12227701	12
Caelo	Paraffin, liquid	12295601	1
Caelo	Paraffin, liquid	12324304	2
GEHE	Paraffin, liquid	12324307	1
Sanacorp	Paraffin, liquid	12324307	1
Caelo	Paraffin, liquid	12324307	22
Phönix Gotha	Paraffin, liquid	12324307	1
Gehe	Paraffin, liquid	12324307	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, liquid	1234307	1
Laborchemie Apolda/Jenn. ...	Paraffin, liquid	12345	1
Caelo	Paraffin, liquid	1235465	1
Fagron	Paraffin, liquid	123586	3
Euro OTC	Paraffin, liquid	125071	1
Anzag	Paraffin, liquid	123586	1
Euro OTC	Paraffin, liquid	125962	1
Caelo	Paraffin, liquid	12H27-N01	2
Fagron	Paraffin, liquid	12H27-N01	8
Fagron	Paraffin, liquid	12H27-N01	1
Fagron	Paraffin, liquid	12J23-N03	4
Caelo	Paraffin, liquid	13005205	1
Caelo	Paraffin, liquid	13005805	17
Lamotte	Paraffin, liquid	13005805	1
Caelo	Paraffin, liquid	13005806	1
Cealo	Paraffin, liquid	1300805	1
Hedinger/Phönix	Paraffin, liquid	1302C008-279571	1
Euro OTC	Paraffin, liquid	1303031-01	10
Noweda	Paraffin, liquid	1303031-01	1
Euro OTC	Paraffin, liquid	130303101	1
Caelo	Paraffin, liquid	1304010-01	1
Euro OTC	Paraffin, liquid	1304010-01	6
Euro OTC	Paraffin, liquid	130401001	1
Caelo	Paraffin, liquid	13050102	22
PHÖNIX 05.06.2013	Paraffin, liquid	13050102	1
Euro OTC	Paraffin, liquid	1307007-01	15
Caelo	Paraffin, liquid	1307007-01	1
Eigenherstellung	Paraffin, liquid	130712HML01	1
Euro OTC	Paraffin, liquid	1309022-01	2
Euro OTC	Paraffin, liquid	130902201	1
Caelo	Paraffin, liquid	1309022-02	2
Spinnrad	Paraffin, liquid	13103308	1
Caelo	Paraffin, liquid	13103309	20
Gehe	Paraffin, liquid	13103309	1
EuRho / Noweda	Paraffin, liquid	131102401	1
Euro OTC	Paraffin, liquid	1311024/01	1
Euro OTC	Paraffin, liquid	1311024-01	13
Caelo	Paraffin, liquid	1311024-01	2
Euro OTC	Paraffin, liquid	131102401	2
Caelo	Paraffin, liquid	13131705	17
Caelo	Paraffin, liquid	13131704	1
Phönix	Paraffin, liquid	13131705	1
Caelo	Paraffin, liquid	1313705	2
Caelo	Paraffin, liquid	1322902	1
Henry Lamotte Oils GmbH...	Paraffin, liquid	13148606	1
Caelo	Paraffin, liquid	13232901	15
Caelo	Paraffin, liquid	13232901ok	1
Caelo	Paraffin, liquid	13232902	14
Caelo	Paraffin, liquid	13232903	8
Caelo	Paraffin, liquid	13248006	2
Caelo	Paraffin, liquid	13248007	8
Euro OTC	Paraffin, liquid	133888	1
Caelo	Paraffin, liquid	13390601	7
Caelo	Paraffin, liquid	13390602	15
Cealo	Paraffin, liquid	13390602	1
Gehe	Paraffin, liquid	13390602	1
Sanacorp	Paraffin, liquid	13390602	1
Devesa/Kehr	Paraffin, liquid	13390602	1
Caelo	Paraffin, liquid	13390604	2
Fagron	Paraffin, liquid	13A3010	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, liquid	13390606	1
Fagron	Paraffin, liquid	13A30-N11	14
Fargon	Paraffin, liquid	13A30-N11	1
Phönix	Paraffin, liquid	13A30-N11	2
Phoenix	Paraffin, liquid	13A30-N11	1
Fagron	Paraffin, liquid	13a30-n11	1
Fagron	Paraffin, liquid	13C19-N09	3
Pönix 10.06.2013	Paraffin, liquid	13C19-N09	1
Großhandel	Paraffin, liquid	13C19-N09	1
Fagron	Paraffin, liquid	13E28-N02	13
Sanacorp	Paraffin, liquid	13E28-N02	1
Fagron	Paraffin, liquid	13G24-N12	6
Fagron	Paraffin, liquid	13G24N12	1
Caelo	Paraffin, liquid	13G24-N12	2
Gehe	Paraffin, liquid	13G24-N12	1
Fagron	Paraffin, liquid	13H02-B02	1
Fagron	Paraffin, liquid	13L17-B01-289491	7
Caelo	Paraffin, liquid	13L17-B01	1
Fagron	Paraffin, liquid	13L17-B01	2
Fagron	Paraffin, liquid	13L17-B01-289492	4
Fiebig	Paraffin, liquid	13L17-B01-289492	1
Fagron	Paraffin, liquid	13L17-B01-289492	1
Fagron	Paraffin, liquid	13L17-b01-289492	1
Fagron	Paraffin, liquid	13L17-B01-289493	4
Fagron	Paraffin, liquid	13L-B01-289491	1
Caelo	Paraffin, liquid	14011303	1
Caelo	Paraffin, liquid	14011514	1
Euro OTC	Paraffin, liquid	1404009-01	14
Caelo	Paraffin, liquid	1404009-01	1
Euro OTC	Paraffin, liquid	140400901	1
Noweda	Paraffin, liquid	1404009-09	1
Euro OTC	Paraffin, liquid	1404010-01	2
Euro OTC	Paraffin, liquid	140401001	2
Fagron	Paraffin, liquid	1404010-01	1
Caelo	Paraffin, liquid	14044504	2
Caelo	Paraffin, liquid	14044507	7
Caelo	Paraffin, liquid	14059701	7
Phönix	Paraffin, liquid	14059701	1
Phoenix 10.10.14	Paraffin, liquid	14059701	1
Caelo	Paraffin, liquid	14059703	4
Caelo	Paraffin, liquid	14059705	1
Caelo	Paraffin, liquid	14059708	4
Caesar&loretz,Anzag	Paraffin, liquid	14059708	1
Euro OTC	Paraffin, liquid	14059708	1
Euro OTC	Paraffin, liquid	1407042-01	13
Caelo	Paraffin, liquid	1407042-01	1
Euro/ Sanacorp	Paraffin, liquid	1409007-01	1
Euro OTC	Paraffin, liquid	140900701	1
Euro OTC	Paraffin, liquid	1412010-01	12
Caelo	Paraffin, liquid	14134601	1
Caelo	Paraffin, liquid	14141406	1
Caelo	Paraffin, liquid	14141408	3
Caelo	Paraffin, liquid	14163804	6
Caelo	Paraffin, liquid	14163807	21
Bombastus	Paraffin, liquid	14187601	1
Caelo	Paraffin, liquid	14202507	9
Caelo	Paraffin, liquid	14202509	12
Caelo	Paraffin, liquid	14305904	6
Caelo	Paraffin, liquid	14305905	8
Caelo	Paraffin, liquid	14357201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, liquid	14357203	9
Noweda 19.05.2015	Paraffin, liquid	14357203	1
Caelo	Paraffin, liquid	14357204	12
Caelo	Paraffin, liquid	14357205	3
Fagron	Paraffin, liquid	14357205	1
Caelo	Paraffin, liquid	1437205	1
Fagron	Paraffin, liquid	14D24-B05	5
Caelo	Paraffin, liquid	14D24-B05	3
Fagron	Paraffin, liquid	14D24-B05-295014	3
Caelo	Paraffin, liquid	14D24-B05-295022	1
Fagron	Paraffin, liquid	14D24-B05-295022	2
Caelo	Paraffin, liquid	14D24B05-295022	1
Fagron	Paraffin, liquid	14d24-b05-295023	1
Hansen & Rosenthal KG/N...	Paraffin, liquid	14D24-B05-295023	1
Fagron	Paraffin, liquid	14D24-B05-295023	3
WE:Phoenix 06.10.2016 ...	Paraffin, liquid	14D24-B05-295023	1
Devesa/ Phönix	Paraffin, liquid	14D24-B05-295022	1
Fagron	Paraffin, liquid	14D24-B05-295023	1
Euro OTC	Paraffin, liquid	1501016-01	5
Euro Pharma	Paraffin, liquid	150101601	1
Euro OTC	Paraffin, liquid	150101601	1
Euro OTC	Paraffin, liquid	1503046-01	8
Eigenherstellung	Paraffin, liquid	150114HML01	1
Caelo	Paraffin, liquid	1503046-01	2
Eu Rho/ Gehe	Paraffin, liquid	1503046-01	1
Eigenherstellung	Paraffin, liquid	150413HML01	1
Euro OTC	Paraffin, liquid	1506005-01	9
Caelo	Paraffin, liquid	150600501	1
Bombastus	Paraffin, liquid	150600501	1
EuRho/Phönix	Paraffin, liquid	1506005-01	1
Euro OTC	Paraffin, liquid	1509011-01	7
anzag	Paraffin, liquid	1509011-01	1
Caelo	Paraffin, liquid	151060104	1
Caelo	Paraffin, liquid	15106101	5
Caelo	Paraffin, liquid	15106102	2
Caelo	Paraffin, liquid	15106104	9
Caelo	Paraffin, liquid	15106107	10
Caelo	Paraffin, liquid	15106109	11
Euro OTC	Paraffin, liquid	1511029-01	4
Euro OTC	Paraffin, liquid	1511030-01	14
Caelo	Paraffin, liquid	15230601	1
Caelo	Paraffin, liquid	15234502	9
Caesar&Loretz GmbH	Paraffin, liquid	15234502	1
Caelo	Paraffin, liquid	15234505	1
Caelo	Paraffin, liquid	15234506	3
Caelo	Paraffin, liquid	15293201	11
Caelo	Paraffin, liquid	15293202	10
Phoenix	Paraffin, liquid	15293202	2
Caelo	Paraffin, liquid	15293203	14
Caelo	Paraffin, liquid	15293204	8
Caelo	Paraffin, liquid	15293205	16
Anzag	Paraffin, liquid	15293205	1
Caesar&Loretz GmbH,AHC,...	Paraffin, liquid	15293205	1
Caelo	Paraffin, liquid	15293209	1
Caelo	Paraffin, liquid	15336702	1
Caelo	Paraffin, liquid	15299201	1
Caelo	Paraffin, liquid	15336703	10
Noweda	Paraffin, liquid	15336703	1
Caelo	Paraffin, liquid	15336704	8
Calo	Paraffin, liquid	15429603	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, liquid	15429603	6
Caelo	Paraffin, liquid	1593205	1
Fagron	Paraffin, liquid	15B16-B03-304611	6
Caelo	Paraffin, liquid	16/080/P	1
Caelo	Paraffin, liquid	1601312050	1
Euro OTC	Paraffin, liquid	16020026-01	1
Euro OTC	Paraffin, liquid	1602026-01	13
Caelo	Paraffin, liquid	1602026-01	1
Euro OTC	Paraffin, liquid	1604009-01	7
Caelo	Paraffin, liquid	1604009-01	1
Bombastus	Paraffin, liquid	1604009-01	1
Euro OTC	Paraffin, liquid	1604045-01	4
Caelo	Paraffin, liquid	1604045-01	1
Euro OTC	Paraffin, liquid	1606032-01	11
Bombastus	Paraffin, liquid	1606032-01	1
Caelo	Paraffin, liquid	16061203	1
Caelo	Paraffin, liquid	1607039-01	1
Euro OTC	Paraffin, liquid	1607039-01	6
Bombastus	Paraffin, liquid	1607039-01	1
Euro OTC	Paraffin, liquid	160703901	1
Euro OTC	Paraffin, liquid	1608007-01	2
Caelo	Paraffin, liquid	160820907	1
Caelo	Paraffin, liquid	160822907	1
Caelo	Paraffin, liquid	160829	1
Caelo	Paraffin, liquid	16082904	12
Caesar&Loretz GmbH,Phön...	Paraffin, liquid	16082904	1
AHD	Paraffin, liquid	16082904	1
Caelo	Paraffin, liquid	16082905	9
Fischar	Paraffin, liquid	16082905	1
Phoenix	Paraffin, liquid	16082906	1
Caelo	Paraffin, liquid	16082906	9
Phoenix	Paraffin, liquid	16082907	1
Caelo	Paraffin, liquid	16082907	8
Caelo	Paraffin, liquid	16082908	8
Euro OTC	Paraffin, liquid	1610017-01	4
EuRho	Paraffin, liquid	1610017-01	1
Phoenix	Paraffin, liquid	1610017-01	1
Caelo	Paraffin, liquid	16149801	12
05.11.2016 Phoenix 2,66EUR	Paraffin, liquid	16149801	1
Caelo	Paraffin, liquid	16149802	19
Ebert KKH 30.09.16	Paraffin, liquid	16149802	1
Caelo	Paraffin, liquid	16149806	2
Caelo	Paraffin, liquid	16171001	6
Caelo	Paraffin, liquid	16171002	1
Caelo	Paraffin, liquid	16238802	6
Euro OTC	Paraffin, liquid	16238802	1
Caelo	Paraffin, liquid	16171008	1
Noweda	Paraffin, liquid	16238802	1
Sanacorp	Paraffin, liquid	16238802	1
Caelo	Paraffin, liquid	16238803	5
Caelo	Paraffin, liquid	16238804	1
Caelo	Paraffin, liquid	16238806	19
Phönix	Paraffin, liquid	16238806	1
Allianz Healthcare	Paraffin, liquid	16238806	1
Caelo	Paraffin, liquid	16238807	3
Caelo	Paraffin, liquid	16238808	4
Caelo	Paraffin, liquid	16238808(1)	1
Caelo	Paraffin, liquid	16238808(2)	1
Caelo	Paraffin, liquid	16238812	1
Caelo	Paraffin, liquid	16238808(3)	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, liquid	16238809	1
Caelo	Paraffin, liquid	16238808(4)	1
Caelo	Paraffin, liquid	16238813	3
Cealo/Gehe	Paraffin, liquid	16238813	1
Caelo	Paraffin, liquid	16238815	1
Caelo	Paraffin, liquid	16238860	1
Hansen & Rosenthal/Noweda	Paraffin, liquid	16A27-B01	1
Fagron	Paraffin, liquid	16A27-B01-323777	3
Caelo	Paraffin, liquid	16A27-B01-323777	1
Hochhaus-Apotheke	Paraffin, liquid	16A27-B01-323777	1
Fagron	Paraffin, liquid	16E30-B05-324041	1
Euro OTC	Paraffin, liquid	1701033-01	4
Euro OTC	Paraffin, liquid	1701052-01	1
Euro/Fiebig	Paraffin, liquid	1704023-01	1
Euro OTC	Paraffin, liquid	1704023-01	1
Euro OTC	Paraffin, liquid	170608P/1610017-01	1
Euro OTC	Paraffin, liquid	17047603	1
Caelo	Paraffin, liquid	17062001	3
Euro OTC	Paraffin, liquid	1704023-0108598729	1
Caelo	Paraffin, liquid	17062002	1
Sanacorp WE:27.12.17 ...	Paraffin, liquid	17062008	1
Euro OTC	Paraffin, liquid	1710025-01	1
Caelo	Paraffin, liquid	17321501	1
Caelo	Paraffin, liquid	17321503	1
Caelo	Paraffin, liquid	1811Q-02050	1
Caelo	Paraffin, liquid	19021302	1
Caelo	Paraffin, liquid	1913Q-02050	1
fargon	Paraffin, liquid	196310	1
Caelo	Paraffin, liquid	2103312050	1
Henry Lamotte Olis	Paraffin, liquid	2x5367601	1
Caelo	Paraffin, liquid	20210312-4	1
Eigenherstellung	Paraffin, liquid	210714H01	1
Euro OTC	Paraffin, liquid	22101402	1
Caelo	Paraffin, liquid	2213A-02050	1
Caelo	Paraffin, liquid	221108111F30-NO	1
Eigenherstellung	Paraffin, liquid	240715HML01	1
Caelo	Paraffin, liquid	25092015	1
Euro OTC	Paraffin, liquid	261113PML01	1
Eigenherstellung	Paraffin, liquid	261113HML01	1
Eigenherstellung	Paraffin, liquid	260115ML01	1
Caelo	Paraffin, liquid	27011405	1
Bombastus	Paraffin, liquid	274958	6
Pönix 10.05.2013	Paraffin, liquid	274958	1
A.H. 03.06.13 ek: 6,15EUR	Paraffin, liquid	274958	1
Bombastus	Paraffin, liquid	281540	3
Bombastus	Paraffin, liquid	284273	3
Bombastus	Paraffin, liquid	288827	2
Bombastus	Paraffin, liquid	290824	2
Eigenherstellung	Paraffin, liquid	291012HJK01	1
Euro OTC	Paraffin, liquid	2911028-01	1
Fagron	Paraffin, liquid	29137	1
Bombastus	Paraffin, liquid	292137	1
Bombastus	Paraffin, liquid	293385	2
Bombastus	Paraffin, liquid	294568	1
Bombastus	Paraffin, liquid	295576	3
Bombastus	Paraffin, liquid	298396	2
Bombastus	Paraffin, liquid	300061	2
Noweda/Bomnastus	Paraffin, liquid	303806	2
Bombastus	Paraffin, liquid	303806	3
Gehe	Paraffin, liquid	310779	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, liquid	3213Q2050	1
fragon	Paraffin, liquid		1
Caelo	Paraffin, liquid	3613A-02050	1
Caelo	Paraffin, liquid	4112M-02048	1
Fagron	Paraffin, liquid	4413M2050	1
Phönix/Henry	Paraffin, liquid	5367601	7
Henry Lamotte Oils GmbH	Paraffin, liquid	5098101	1
Henry Lamotte/Anzag	Paraffin, liquid	5367601	3
Henry Lamotte	Paraffin, liquid	5367601	14
Henry Lamotte GmbH, Mer...	Paraffin, liquid	5367601	2
Phönix/Lamotte	Paraffin, liquid	5367601	2
Henry Lamotte/GEHE	Paraffin, liquid	5367601	2
Kehr	Paraffin, liquid	5367601	1
Caelo	Paraffin, liquid	5367601	2
Alliance Healthcare	Paraffin, liquid	5367601	2
Lamotte	Paraffin, liquid	5367601	4
Lamotte/GEHE	Paraffin, liquid	5367601	1
Henry Lamotte Olis GmbH	Paraffin, liquid	5367601	2
Sanacorp	Paraffin, liquid	5367601	1
LAMOTTE/NOWEDA	Paraffin, liquid	5367601	1
Henry Lamotte Oils/Phönix	Paraffin, liquid	5367601	1
Henry Lamotte Oils GmbH	Paraffin, liquid	5367601	1
H.Lamotte/Noweda	Paraffin, liquid	5367601	1
Hedinger	Paraffin, liquid	5367601	1
Henry Lamotte	Paraffin, liquid	5593601	8
Noweda / Henry Lamotte	Paraffin, liquid	5593601	4
Henry Lamotte Olis GmbH...	Paraffin, liquid	5593601	2
lamotte noweda,	Paraffin, liquid	5593601	1
Caelo	Paraffin, liquid	5593601	1
Lamotte/Gehe	Paraffin, liquid	5593601	1
Henry lamotte Oils GmbH...	Paraffin, liquid	5593601	1
Noweda / Henry Lamotte ...	Paraffin, liquid	5593601	1
Sanacorp	Paraffin, liquid	5593601	2
Henry Lamotte Oils GmbH	Paraffin, liquid	5593601	5
Henry Lamotte oils GmbH...	Paraffin, liquid	5593601	5
Lamotte / Noweda	Paraffin, liquid	5593601	1
Noweda/Lamotte	Paraffin, liquid	5593601	1
Alliance health	Paraffin, liquid	5593601	1
Henry Lamotte/Sana	Paraffin, liquid	5593601	2
Lamotte / Alliance Heal...	Paraffin, liquid	5593601	1
Caelo	Paraffin, liquid	57/AD1	1
Caelo	Paraffin, liquid	58/AKD55	1
Henry Lamotte	Paraffin, liquid	5863401	3
henry lamotte oils/Sana...	Paraffin, liquid	5863401	1
Caelo	Paraffin, liquid	59/AKD2	1
Sanacorp	Paraffin, liquid	59/AKD4	1
Henry Lamotte/Phönix	Paraffin, liquid	7316019	1
Caelo	Paraffin, liquid	83054508	1
Henry Lamotte GmbH, Mer...	Paraffin, liquid	79460	1
Fagron	Paraffin, liquid	959561	4
Phoenix	Paraffin, liquid	959561	1
Fagron	Paraffin, liquid	B05-295022	1
Euro OTC	Paraffin, liquid	Ch1211026-01	1
Eigenherstellung	Paraffin, liquid		1
ACM	Paraffin, liquid	KNC2106J14	1
ACM	Paraffin, liquid	KNC2160J14	1
Euro OTC	Paraffin, liquid	1102009-01	1
Herbasinica	Paraffin, liquid	11354102	1
Fagron	Paraffin, liquid	1198738	1
Sanacorp	Paraffin, liquid	1211027-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
EUro/AHD	Paraffin, liquid	1211027-01	1
Caelo	Paraffin, liquid	12324306	1
Sanacorp	Paraffin, liquid	1304010-01	1
Celoa vdl	Paraffin, liquid	131102401	1
Gutenberg-Apotheke	Paraffin, liquid	13248007	1
Henry Lamotte/Sanacorp	Paraffin, liquid	1313703	1
Fagron	Paraffin, liquid	13E28-N02	1
Caeo	Paraffin, liquid	14305905	1
Euro OTC	Paraffin, liquid	1506021-01	1
Euro	Paraffin, liquid	1511030-01	1
HansenRosenthalKG/Hansen	Hydrophobic base gel DAC	1078667	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes®	13A089	1
Ichthyol	Unguentum Cordes®	13A091	1
Caelo	Unguentum Cordes®	13A113	1
Ichthyol Gesellschaft	Unguentum Cordes®	13A171	1
Ichthyol-Gesellschaft	Unguentum Cordes®	13A172	1
Ichthyol Gesellschaft/AHD	Unguentum Cordes®	13A180	1
Fagron	Wool alcohol ointment DAB	13D09-T04-007315	1
Ichthyol/gehe	Unguentum Cordes®	14A005	1
Ichthyol/gehe	Unguentum Cordes®	14A099	1
Caelo	Unguentum Cordes®	14A095	1
Ichthyol/Phönix	Unguentum Cordes®	15A020	2
AHD	Unguentum Cordes®	15A021	1
Ichthyol/Phönix	Unguentum Cordes®	15A023	1
A.H.D.	Unguentum Cordes®	15A060	1
Ichthyol/Phönix	Unguentum Cordes®	15A029	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	15A042	1
Ichthyol/Noweda	Unguentum Cordes®	15A060	1
Ichthyol/Phönix	Unguentum Cordes®	15A062	1
Ichthyol/ AHD	Unguentum Cordes®	15A091	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	15A091	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes®	15A123	1
Caelo	Unguentum Cordes®	15A123	1
Ichthyol/Phönix	Unguentum Cordes®	15A124	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes®	15A124	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	15A156	1
Ichthyol/Phönix	Unguentum Cordes®	15A162	1
Ichthyol/Gehe	Unguentum Cordes®	15A78	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes®	15A191	1
Caelo	Wool alcohol ointment DAB	16331702	8
Ichthyol/Phönix	Unguentum Cordes®	16A023	1
Cordes	Unguentum Cordes®	16A046	1
Noweda	Unguentum Cordes®	16A047	1
Sanacorp, 15.08.2016, 3...	Unguentum Cordes®	16A048	1
ICHTHYOL-GESELLSCHAFT C...	Unguentum Cordes®	16A083	1
Galderma/Phönix	Unguentum Cordes®	16A077	1
Ichthyol / Noweda	Unguentum Cordes®	16A123	1
ICHTHYOL-GESELLSCHAFT/P...	Unguentum Cordes®	16A173	1
ICHTHYOL-GESELLSCHAFT/P...	Unguentum Cordes®	16A140	1
Ichthyol-Gesellschaft	Unguentum Cordes®	17A126	1
Ichthyol-Gesellschaft	Unguentum Cordes®	17A003	1
Gehe	Unguentum Cordes®	17A004	1
Ichthyol	Unguentum Cordes®	17A004	1
Ichthyol-Gesellschaft /...	Unguentum Cordes®	17A004	2
Cordes	Unguentum Cordes®	17A004	1
Ichthyol/Noweda	Unguentum Cordes®	17A004	1
Caelo	Unguentum Cordes®	17A004	1
ICHTYOL-Gesellsch./Sana...	Unguentum Cordes®	17A004	1
Ichthyol/Noweda	Unguentum Cordes®	17A005	1
Ichthyol Gesells./AHD	Unguentum Cordes®	17A005	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/Noweda	Unguentum Cordes [®]	17A007	1
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	17A007	1
Fagron	Unguentum Cordes [®]	17A005	1
Ichthyol Gesellschaft/AHD	Unguentum Cordes [®]	17A054	1
ICHTHYOL-GESELLSCHAFT	Unguentum Cordes [®]	17A049	1
Ichthyol-Gesellschaft	Unguentum Cordes [®]	17A049	1
Noweda-Ichthyol	Unguentum Cordes [®]	17A054	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A056	1
Ichthyol	Unguentum Cordes [®]	17A056	1
Ichthyol/Gehe	Unguentum Cordes [®]	17A057	1
Ichthyol	Unguentum Cordes [®]	17A055	1
Ichthyol/ Phoenix	Unguentum Cordes [®]	17A069	1
Ichthyol / Phönix	Unguentum Cordes [®]	17A098	1
Ichthyol	Unguentum Cordes [®]	17A098	1
Phönix W.E:09.03.18 E...	Unguentum Cordes [®]	17A098	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A104	1
Phoenix	Unguentum Cordes [®]	17A104	1
Ichthyol-Noweda	Unguentum Cordes [®]	17A105	1
ICHTHYOL / 20	Unguentum Cordes [®]	17A106	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	17A106	1
Cordes/Noweda	Unguentum Cordes [®]	17A106	1
Ichthyol/Phönix	Unguentum Cordes [®]	17A112	2
Apomix	Eye ointment, emulsifying DAC	17A111113	1
Ichthyol	Unguentum Cordes [®]	17A113	2
Ichthyol-Gesellschaft C...	Unguentum Cordes [®]	17A113	1
Phoenix	Unguentum Cordes [®]	17A112	1
Ichthyol/ Phoenix	Unguentum Cordes [®]	17A126	1
Ichthyol/GEHE	Unguentum Cordes [®]	17A127	1
Apomix	Eye ointment, emulsifying DAC	18B021214	1
Fagron	Wool alcohol ointment DAB	23033310	1
Bombastus	Wool alcohol ointment DAB	279928	4
Bombastus	Wool alcohol ointment DAB	283619	3
Bombastus	Wool alcohol ointment SR wit...	285952	2
Bombastus	Hydrophilic ointment DAB	286543	1
Bombastus	Wool alcohol ointment DAB	287619	2
Bombastus	Wool alcohol ointment DAB	290735	2
Bombastus	Wool alcohol ointment SR wit...	291177	2
Bombastus	Wool alcohol ointment DAB	291795	5
Beiersdorf/Gehe	Wool alcohol ointment DAB	40726610	1
Beiersdorf	Wool alcohol ointment DAB	42541910	1
Lamotte/Anzag	Hydrophilic ointment DAB	5112A-02539	1
Lamotte	Hydrophilic ointment DAB	5651600	2
HENRY LAMOTTE OILS GmbH	Hydrophilic ointment DAB	5860400	2
Beiersdorf	Wool alcohol ointment DAB	11235310	1
Caelo	Wool alcohol ointment DAB	13242603	6
Beiersdorf/AHD	Wool alcohol ointment DAB	22412510	1
Beiersdorf	Wool alcohol ointment DAB	24655410	1
Gehe	Wool alcohol ointment SR wit...	276135	1
Bombastus	Wool alcohol ointment DAB	277022	1
Bombastus	Wool alcohol ointment SR wit...	277022	2
Sanacorp	Wool alcohol ointment DAB	277190	2
Bombastus	Wool alcohol ointment DAB	279105	3
Sanacorp	Wool alcohol ointment DAB	279927	1
Bombastus	Wool alcohol ointment SR wit...	280058	2
Bombastus	Wool alcohol ointment SR wit...	280059	1
Bombastus	Wool alcohol ointment SR wit...	283176	4
Beiersdorf, Phönix	Wool alcohol ointment DAB	30721310	1
Beiersdorf, Phönix	Wool alcohol ointment DAB	31121710	1
Beiersdorf	Wool alcohol ointment DAB	31121710	1
Beiersdorf	Wool alcohol ointment DAB	31241810	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Ichthyol/Anzag	Unguentum Cordes [®]	12A055	1
Sanova	Unguentum Cordes [®]	11A040	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	12A006	1
Gehe	Unguentum Cordes [®]	11A161	1
Ichthyol/N/14.02.13	Unguentum Cordes [®]	12A131	1
Caelo	Unguentum Cordes [®]	12A140	1
ichthyol	Unguentum Cordes [®]	12A115	1
Caelo	Unguentum Cordes [®]	12A141	1
Ichtyol-Gesellschaft	Unguentum Cordes [®]	12A184	1
Ichthyol Jenne	Unguentum Cordes [®]	12A211	1
Caelo	Unguentum Cordes [®]	13A004	1
Ichthyol/Phönix	Unguentum Cordes [®]	13A010	1
Sanacorp	Unguentum Cordes [®]	13A042	1
Ichthyol/Sanacorp	Unguentum Cordes [®]	13A010	1
ichthyol/N/PZN:01062983	Unguentum Cordes [®]	13A050	1
Ichthyl Phönix	Unguentum Cordes [®]	13A053	1
Ichthyol	Unguentum Cordes [®]	13A109	1
Ichthyol/gehe	Unguentum Cordes [®]	13A114	1
ichthyol/AHD	Unguentum Cordes [®]	13A109	1
Ichthyol/GEHE	Unguentum Cordes [®]	13a138	1
Bombastus	Hydrophilic ointment DAB	269534	2
Anzag	Hydrophilic ointment DAB	271047	1
Bombastus	Hydrophilic ointment DAB	274003	1
Bombastus	Hydrophilic ointment DAB	274002	1
Caelo	Hydrophilic ointment DAB	13171102	1
Bombastus	Hydrophilic ointment DAB	276670	1
Bombastus	Hydrophilic ointment DAB	279000	2
Apomix	Eye ointment, emulsifying DAC	5B220413	1
Apomix	Eye ointment, emulsifying DAC	15A241013	1
Noweda	Paraffin, white	255510	1
Bombastus	Paraffin, white	258194	1
Bombastus	Paraffin, white	258551	1
Gehe	Paraffin, white	261731	1
Bombastus	Paraffin, white	263683	1
Bombastus	Paraffin, white	263759	1
Noweda	Paraffin, white	263762	1
Bombastus	Paraffin, white	265508	1
Bombastus	Paraffin, white	265814	2
Bombastus	Paraffin, white	267010	2
Bombastus	Paraffin, white	268016	1
Bombastus	Paraffin, white	268099	2
Gehe	Paraffin, white	268099	1
Bombastus	Paraffin, white	268390	1
Gehe	Paraffin, white	268391	1
Bombastus	Paraffin, white	269731	2
Bombastus	Paraffin, white	269732	3
Bombastus	Paraffin, white	269733	3
Bombastus	Paraffin, white	269734	5
Bombastus	Paraffin, white	270163	2
Noweda	Paraffin, white	270163	3
Bombastus	Paraffin, white	270165	2
Bombastus	Paraffin, white	270489	3
Bombastus	Paraffin, white	270490	1
Bombastus	Paraffin, white	270491	1
Phoenix	Paraffin, white	270491	1
Anzag 08.10.12	Paraffin, white	270491	1
Bombastus	Paraffin, white	270633	1
Phönix Gotha	Paraffin, white	270633	1
Bombastus	Paraffin, white	270634	2
Phönix Gotha	Paraffin, white	270634	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Paraffin, white	270635	2
Bombastus	Paraffin, white	270834	4
Bombastus	Paraffin, white	270976	5
Bombastus	Paraffin, white	270835	1
Bombastus	Paraffin, white	270977	1
Bombastus	Paraffin, white	270978	3
Bombastus	Paraffin, white	270979	6
Bombastus	Paraffin, white	270982	3
Bombastus	Paraffin, white	270983	1
Bombastus	Paraffin, white	271213	4
Bombastus	Paraffin, white	271260	1
Bombastus	Paraffin, white	271261	2
Noweda	Paraffin, white	271261	1
Gehe	Paraffin, white	271435	1
Bombastus	Paraffin, white	271435	2
Bombastus	Paraffin, white	271437	5
Anzag	Paraffin, white	271435	1
Caelo	Paraffin, white	1061301	1
Caelo	Paraffin, white	1021302	1
Caelo	Paraffin, white	1081409	1
Caelo	Paraffin, white	2061201	1
Caelo	Paraffin, white	213E-02549	1
Caelo	Paraffin, white	402312550	1
Caelo	Paraffin, white	4031406	1
Caelo	Paraffin, white	3071202	1
Caelo	Paraffin, white	3091209	1
Hedinger/Noweda	Paraffin, white	4-219072	1
Caelo	Paraffin, white	5041303	1
Caelo	Paraffin, white	6061302	1
Caelo	Paraffin, white	5061412	1
Caelo	Paraffin, white	613Q-02549	1
Bombastus	Paraffin, white	6861393	1
Caelo	Paraffin, white	713I-02549	1
Caelo	Paraffin, white	8012014D	1
Caelo	Paraffin, white	8021401	1
Caelo	Paraffin, white	8031302	1
Caelo	Paraffin, white	8051303	1
Caelo	Paraffin, white	8101302	1
Caelo	Paraffin, white	813E-02549	1
Caelo	Paraffin, white	9274706	1
Caelo	Paraffin, white	100/AKD8	1
Caelo	Paraffin, white	10051402	1
Caelo	Paraffin, white	10071304	1
Bombastus	Paraffin, white	1001312549	1
Caelo	Paraffin, white	10071305	1
Caelo	Paraffin, white	10091301	1
Caelo	Paraffin, white	10094704	1
Caelo	Paraffin, white	10155402	1
Caelo	Paraffin, white	10220214	1
Sanacorp	Paraffin, white	102/AKD1	1
Caelo	Paraffin, white	10271001	1
Caelo	Paraffin, white	1101559	1
Caelo	Paraffin, white	11011309	1
Fagron	Paraffin, white	10H12-T01	1
Caelo	Paraffin, white	11040603	1
Caelo	Paraffin, white	11041302	1
Caelo	Paraffin, white	11031403	1
Anzag	Paraffin, white	1105029-01	1
Caelo	Paraffin, white	11060209	3
Caelo	Paraffin, white	11060210	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	11060211	1
Caelo	Paraffin, white	11060212	2
Caelo	Paraffin, white	11060213	6
Euro OTC	Paraffin, white	1108046-01	4
Hansen&Rosenthal/Phönix	Paraffin, white	1108046-01	1
Euro OTC	Paraffin, white	1109002-01	2
Caelo	Paraffin, white	11091402	1
Caelo	Paraffin, white	11104901	1
Caelo	Paraffin, white	11104903	1
Caelo	Paraffin, white	11104904	2
Caelo	Paraffin, white	11104910	1
Caelo	Paraffin, white	11108711	1
Caelo	Paraffin, white	11108710	1
Caelo	Paraffin, white	11190913	1
Caelo	Paraffin, white	11115501	1
Caelo	Paraffin, white	11265701	1
Caelo	Paraffin, white	11265708	2
Caelo	Paraffin, white	11265703	1
vdl	Paraffin, white	11265710	1
Caelo	Paraffin, white	11265711	1
Caelo	Paraffin, white	11283902	1
Caelo	Paraffin, white	11301301	5
Caesar & Loretz GmbH, 4. .	Paraffin, white	11301301	1
Noweda	Paraffin, white	11301301	1
Caelo	Paraffin, white	11315401	2
Caelo	Paraffin, white	11315402	1
Gehe	Paraffin, white	11315403	1
Caelo	Paraffin, white	11315408	3
Caelo	Paraffin, white	11315412	1
Caelo	Paraffin, white	11315413	2
Caelo	Paraffin, white	11374901	1
Caelo	Paraffin, white	11374903	4
Caelo	Paraffin, white	11374906	6
Noweda	Paraffin, white	11374906	1
Fagron	Paraffin, white	11G11-T02	2
Fagron	Paraffin, white	11K14-T01	7
Fagron	Paraffin, white	12/17-T01	2
Phönix 08.04.2013	Paraffin, white	12/17-T01	1
Fagron	Paraffin, white	12/17-To1	1
Caesar&loretz,Anzag	Paraffin, white	12003601	2
Caelo	Paraffin, white	12003603	1
Euro OTC	Paraffin, white	1201035-01	1
Caelo	Paraffin, white	12012604	5
Caelo	Paraffin, white	12012605	1
Caelo	Paraffin, white	12012609	8
Caelo	Paraffin, white	12012707	2
Sanacorp	Paraffin, white	12031506	1
Gehe	Paraffin, white	12032404	1
Caelo	Paraffin, white	12031202	1
Caelo	Paraffin, white	12032404	3
Caelo	Paraffin, white	12032405	1
Caelo	Paraffin, white	12032412	2
Erika apotheke	Paraffin, white	12032413	1
Caelo	Paraffin, white	12032413	21
ANZAG	Paraffin, white	12032413	1
Phönix Gotha	Paraffin, white	12032413	1
Phönix	Paraffin, white	12032413	1
Caelo	Paraffin, white	12033809	1
Caelo	Paraffin, white	1203412	1
Caelo	Paraffin, white	1203413	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Paraffin, white	1205048-01	1
Euro OTC	Paraffin, white	1205048-01	3
Caelo	Paraffin, white	12071804	1
Caelo	Paraffin, white	12061201	1
Caelo	Paraffin, white	12071805	5
Caesar & Loretz GmbH, 4...	Paraffin, white	12071805	1
Caelo	Paraffin, white	12081404	1
Caelo	Paraffin, white	12078105	1
Caelo	Paraffin, white	12081601	3
Euro OTC	Paraffin, white	1210006-01	9
Euro OTC	Paraffin, white	1210009-01	3
Euro OTC	Paraffin, white	1210010-01	2
Caelo	Paraffin, white	12111203	1
Caelo	Paraffin, white	12111307	1
Großhandel	Paraffin, white	12115101	1
Caelo	Paraffin, white	12115101	10
Anzag 27.03.13 ek: 3,66EUR	Paraffin, white	12115101	1
Caelo	Paraffin, white	12115102	17
Caelo	Paraffin, white	12115105	1
Caelo	Paraffin, white	12115106	7
Noweda	Paraffin, white	12115106	1
Gehe	Paraffin, white	12115106	1
GEHE	Paraffin, white	12115106	1
Doroderm GmbH/ Noweda	Paraffin, white	12115106	2
Sanacorp	Paraffin, white	12115107	1
Caelo	Paraffin, white	12115107	10
Fiebig 22.02.13 ek:3,63EUR	Paraffin, white	12115107	1
Phönix	Paraffin, white	12115107	1
Caelo	Paraffin, white	12115108	6
Fagron	Paraffin, white	12115107	1
Caelo	Paraffin, white	12115113	1
Caelo	Paraffin, white	12121201	1
Caelo	Paraffin, white	12121204	1
Sanacorp	Paraffin, white	1215107	1
Caelo	Paraffin, white	12182504	1
Caelo	Paraffin, white	12182505	9
Anzag/	Paraffin, white	12182505	2
Noweda	Paraffin, white	12182505	1
Caelo	Paraffin, white	12190401	4
Sanacorp	Paraffin, white	12190401	1
Caelo	Paraffin, white	12209203	3
Caelo	Paraffin, white	122092004	1
Noweda	Paraffin, white	12209203	1
Großhandel	Paraffin, white	12209204	1
Caelo	Paraffin, white	12209204	18
Caesar & Loretz GmbH, 4...	Paraffin, white	12209204	1
Anzag 27.12.12 ek:5.45	Paraffin, white	12209204	1
Anzag	Paraffin, white	12209204	1
Caelo	Paraffin, white	12209210	29
Noweda	Paraffin, white	12209210	1
Anzag 15.03.13 ek: 5,61EUR	Paraffin, white	12209210	1
Phönix	Paraffin, white	12209210	1
Caesar & Loretz GmbH, 4...	Paraffin, white	12209210	1
Pönix 19.04.2013	Paraffin, white	12209210	2
Cealo	Paraffin, white	12209210	1
Caelo	Paraffin, white	12306202	2
Caelo	Paraffin, white	12306302	6
Alliance Healthcare	Paraffin, white	12306302	1
Caelo	Paraffin, white	12306303	1
Caelo	Paraffin, white	12306401	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	1231010	1
Caelo	Paraffin, white	12311003	5
Caelo	Paraffin, white	12311004	5
Gehe	Paraffin, white	12311004	1
Caelo	Paraffin, white	12311010	14
Sanacorp	Paraffin, white	12311010	1
Caelo	Paraffin, white	12311011	10
Caelo	Paraffin, white	12321902	1
Caelo	Paraffin, white	12321905	4
Caesar/Lorenz	Paraffin, white	12321906	1
Caelo	Paraffin, white	12321906	4
Caelo	Paraffin, white	12321908	4
Caelo	Paraffin, white	12337701	8
Caelo	Paraffin, white	12337703	18
Caelo	Paraffin, white	12337704	14
Sanacorp	Paraffin, white	12337704	2
Caesar/Lorenz	Paraffin, white	12337704	1
Caelo	Paraffin, white	12337706	14
Caesar & Loretz GmbH, 4...	Paraffin, white	12337706	1
Gehe am 27.2.14	Paraffin, white	12337706	1
Caelo	Paraffin, white	12337707	10
Euro OTC	Paraffin, white	12337707	1
Noweda	Paraffin, white	12337707	1
Caelo	Paraffin, white	12337708	18
Sanacorp	Paraffin, white	12337708	1
Caesar&Loretz GmbH/Anzag	Paraffin, white	12337708	1
Caelo	Paraffin, white	12337709	16
Phoenix 18.03.2014	Paraffin, white	12337709	1
Caelo	Paraffin, white	12374501	5
Großhandel	Paraffin, white	12374502	1
Caelo	Paraffin, white	12374502	13
AHC 08.07.13 EK5,61	Paraffin, white	12374502	1
Caelo	Paraffin, white	1237706	1
Fagron	Paraffin, white	12B20-T02	1
Phönix 16.10.12	Paraffin, white	12D16-T01	1
Fagron	Paraffin, white	12D16-T01	2
Fagron	Paraffin, white	12D16-To1	1
Fagron	Paraffin, white	12D19-T54-000715	4
Fagron	Paraffin, white	12H10-T05	2
Fagron	Paraffin, white	12I17-T01	2
Fagron	Paraffin, white	12J31-T52-000963	1
Fagron	Paraffin, white	12K01-T52-000717	1
Fagron	Paraffin, white	12K01-T52	1
Fagron	Paraffin, white	12K01-T52-002659	2
Caelo	Paraffin, white	12K01-T52-002659	1
Euro OTC	Paraffin, white	12K01-T52-002659	1
Fagron	Paraffin, white	13/19-T05-006192	1
Caelo	Paraffin, white	13005701	9
Fagron	Paraffin, white	13/19-T05-006193	1
Caelo	Paraffin, white	13005708	2
Caelo	Paraffin, white	13005710	9
Caelo	Paraffin, white	13005711	12
Caelo	Paraffin, white	13005712	13
Caelo	Paraffin, white	13005713	4
Caelo	Paraffin, white	13005719	9
Caelo	Paraffin, white	13005720	2
Caelo	Paraffin, white	13011407	1
Euro OTC	Paraffin, white	1302022-01	8
Caelo	Paraffin, white	13028302	1
Caelo	Paraffin, white	13031301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Paraffin, white	1305031-01	13
Caelo	Paraffin, white	1305031-01	1
Caelo	Paraffin, white	13051305	1
Caelo	Paraffin, white	130510	1
Caelo	Paraffin, white	13084902	14
Caelo	Paraffin, white	13084905	13
Caelo	Paraffin, white	13084903	1
Caesar/Noweda	Paraffin, white	13084905	1
Caelo	Paraffin, white	13084906	14
Sanacorp	Paraffin, white	13084906	3
Phoenix	Paraffin, white	13084906	1
Caelo	Paraffin, white	13084907	1
Caelo	Paraffin, white	13095601	10
Caelo	Paraffin, white	13095602	1
canea, gehe	Paraffin, white	13095701	2
Caelo	Paraffin, white	13095701	4
Caesar & Loretz GmbH Hi...	Paraffin, white	13095701	1
Caelo	Paraffin, white	13109010	1
Caelo	Paraffin, white	1310909	1
Caelo	Paraffin, white	1310911	1
Caelo	Paraffin, white	13121303	1
gehe 215.08.2013	Paraffin, white	1313009	1
Caelo	Paraffin, white	13130901	6
Caelo	Paraffin, white	13130907	6
Caelo	Paraffin, white	13130909	15
Sanacorp	Paraffin, white	13130909	2
gehe 30.09.13	Paraffin, white	13130909	1
Caelo	Paraffin, white	13130910	7
Caesar & Loretz GmbH, 4...	Paraffin, white	13130911	1
Caelo	Paraffin, white	13130911	8
Großhandel	Paraffin, white	13130911	1
Caelo	Paraffin, white	13130912	13
Caelo	Paraffin, white	13163701	6
Caelo	Paraffin, white	13163602	2
Anzag	Paraffin, white	13163701	1
Caelo	Paraffin, white	13163702	4
Caelo	Paraffin, white	1318504	1
Caelo	Paraffin, white	13185502	9
Caelo	Paraffin, white	13185501	1
Caelo	Paraffin, white	13174901	1
Caelo	Paraffin, white	13185504	5
Caelo	Paraffin, white	13185505	1
Caelo	Paraffin, white	13185506	5
Caolo/Anzag	Paraffin, white	13185505	1
Caelo	Paraffin, white	13185507	6
Caelo	Paraffin, white	13185512	9
Caelo	Paraffin, white	13185513	1
Caelo	Paraffin, white	13185514	2
Caelo	Paraffin, white	1326408	1
Caelo	Paraffin, white	13268408	22
Caelo	Paraffin, white	13268409	14
Caelo	Paraffin, white	13268410	9
AHD	Paraffin, white	13268410	2
Phönix	Paraffin, white	13268410	1
Caelo	Paraffin, white	13268411	6
Caelo	Paraffin, white	13268420	4
Sanacorp	Paraffin, white	13268420	1
Großhandel	Paraffin, white	13268421	1
Caelo	Paraffin, white	13268421	12
Noweda	Paraffin, white	13268421	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	13268423	3
Phoenix	Paraffin, white	13268423	1
Caelo	Paraffin, white	13336702	2
Caelo	Paraffin, white	13333101	1
Caelo	Paraffin, white	13338702	5
Noweda	Paraffin, white	13337706	1
Caelo	Paraffin, white	13338703	13
Fagron	Paraffin, white	13338703	1
Phoenix am 10.03.2014	Paraffin, white	13338703	1
Caelo	Paraffin, white	13338704	3
Caelo	Paraffin, white	13338707	7
Caelo	Paraffin, white	13338710	13
Berkel AHK	Paraffin, white	13338710	1
Caelo	Paraffin, white	13338711	1
Sanacorp	Paraffin, white	13338712	1
Caelo	Paraffin, white	13338712	8
Herbasinica	Paraffin, white	13338712	1
Caelo	Paraffin, white	13338713	6
Caelo	Paraffin, white	13338714	7
Caelo	Paraffin, white	13341602	4
Phonix	Paraffin, white	13341604	1
Caelo	Paraffin, white	13341605	7
Caelo	Paraffin, white	13341606	6
	Paraffin, white	13341606	1
Caelo	Paraffin, white	13341608	8
Caelo	Paraffin, white	13341612	8
Caelo	Paraffin, white	13365103	2
Sanacorp 05.07.14 EK: 6...	Paraffin, white	13365104	2
Caelo	Paraffin, white	13365104	3
Caelo	Paraffin, white	13380702	4
Caelo	Paraffin, white	13380703	12
Caelo	Paraffin, white	13380707	11
Caelo	Paraffin, white	13380719	7
Caelo	Paraffin, white	13380708	10
Caelo	Paraffin, white	13380720	24
Caelo	Paraffin, white	13380721	15
Caelo	Paraffin, white	1338714	1
Fagron	Paraffin, white	13D02-T02-002660	1
Caelo	Paraffin, white	13D02-T02-002660	1
Fagron	Paraffin, white	13D02-T02-004621	1
Fagron	Paraffin, white	13I19-T05-006192	2
Fagron	Paraffin, white	13I19-T05-006193	1
Fagron	Paraffin, white	13I19-T05-008452	2
Beiersdorf/Gehe	Paraffin, white	1401017-01	1
Euro OTC	Paraffin, white	1401017-01	3
Euro OTC	Paraffin, white	1401018-01	2
Fagron	Paraffin, white	1401453	1
Anzag	Paraffin, white	1402501	1
Caelo	Paraffin, white	14025101	4
Caelo	Paraffin, white	14025201	12
Euro OTC	Paraffin, white	1405014-01	9
Bombastus	Paraffin, white	1406212549	1
Caelo	Paraffin, white	14081205	1
Caelo	Paraffin, white	14087402	4
Caelo	Paraffin, white	14087403	2
Berkel AHK	Paraffin, white	14087403	1
Caelo	Paraffin, white	14087404	11
Caelo	Paraffin, white	14087410	6
Caelo	Paraffin, white	14087405	3
Caelo	Paraffin, white	14124702	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	14124703	3
Phönix	Paraffin, white	14124703	1
Caelo	Paraffin, white	14124704	6
Caelo	Paraffin, white	14124705	1
Caelo	Paraffin, white	14124706	3
Caelo	Paraffin, white	14124707	14
Caelo	Paraffin, white	14189004	4
Caelo	Paraffin, white	14189002	7
Phönix/Lihua	Paraffin, white	14189006	3
Caelo	Paraffin, white	14189018	7
Caelo	Paraffin, white	14189006	3
Fagron	Paraffin, white	14189019	1
Caelo	Paraffin, white	14200901	5
Fagron	Paraffin, white	14201001	1
Caelo	Paraffin, white	14201001	2
Caelo	Paraffin, white	14301804	7
Caelo	Paraffin, white	14347302	11
Caelo	Paraffin, white	14324301	7
Caelo	Paraffin, white	14355806	10
Caelo	Paraffin, white	14324402	1
Fagron	Paraffin, white	14C14T01	3
Fagron	Paraffin, white	14e30-t07-012095	1
Fagron	Paraffin, white	14J03-T03-020827	5
Fagron	Paraffin, white	14E30-T07-012095	1
Fagron	Paraffin, white	14j03-to3-020827	1
Caelo	Paraffin, white	15031402	1
Caelo	Paraffin, white	15041305	1
Caelo	Paraffin, white	15050002	4
Caelo	Paraffin, white	15050003	1
Caelo	Paraffin, white	15060003	12
Caelo	Paraffin, white	15060005	8
Caelo	Paraffin, white	15060012	4
Caelo	Paraffin, white	15060008	3
Caelo	Paraffin, white	15060013	3
Caelo	Paraffin, white	15060014	12
Caelo	Paraffin, white	15081406	1
Caelo	Paraffin, white	15101308	1
Caelo	Paraffin, white	15162503	7
Caelo	Paraffin, white	15162504	7
Caelo	Paraffin, white	15162505	10
Caelo	Paraffin, white	15162507	8
Caelo	Paraffin, white	15162508	2
Caelo	Paraffin, white	15162512	3
Caelo	Paraffin, white	15204301	9
Caelo	Paraffin, white	15204302	8
Caelo	Paraffin, white	15204303	17
Caelo	Paraffin, white	15204304	6
Caelo	Paraffin, white	15219903	3
Caelo	Paraffin, white	15278801	1
Caelo	Paraffin, white	15278804	6
Caelo	Paraffin, white	15278805	3
Caelo	Paraffin, white	15330401	1
Caelo	Paraffin, white	15330502	1
Caelo	Paraffin, white	15339806	12
Caelo	Paraffin, white	15336810	4
Caelo	Paraffin, white	15336809	6
Caelo	Paraffin, white	15349701	5
Caelo	Paraffin, white	15349708	5
Caelo	Paraffin, white	1544150	1
Caelo	Paraffin, white	15441501	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	15441502	2
Caelo	Paraffin, white	15441510	1
Caelo	Paraffin, white	16029211	9
Caelo	Paraffin, white	16029221	6
Caelo	Paraffin, white	16071306	1
Caelo	Paraffin, white	16029213	1
Caelo	Paraffin, white	16061403	1
Caelo	Paraffin, white	16138003	12
Sanacorp	Paraffin, white	16AKD2	1
Caelo	Paraffin, white	17011309	1
Caelo	Paraffin, white	17081202	1
Caelo	Paraffin, white	19071401	1
Caelo	Paraffin, white	17121206	1
Caelo	Paraffin, white	18041315	1
Caelo	Paraffin, white	19091308	1
Bombastus	Paraffin, white	2002312549	1
Caelo	Paraffin, white	20031303	1
Phönix Gotha	Paraffin, white	200567	1
Caelo	Paraffin, white	20061206	1
Caelo	Paraffin, white	20081304	1
Phönix Gotha	Paraffin, white	201387	1
Bombastus	Paraffin, white	2014Q2549	1
Caelo	Paraffin, white	20160812-4	1
Caelo	Paraffin, white	21101306	1
Caelo	Paraffin, white	22011401	1
Caelo	Paraffin, white	22081206	1
Caelo	Paraffin, white	22111306	1
Caelo	Paraffin, white	23051101	1
Caelo	Paraffin, white	23051404	1
Caelo	Paraffin, white	23071202	1
Caelo	Paraffin, white	25021403	1
Caelo	Paraffin, white	250414002	1
Caelo	Paraffin, white	26061307	1
Caelo	Paraffin, white	26081409	1
Caelo	Paraffin, white	26091303	1
Caelo	Paraffin, white	26101202	1
Caelo	Paraffin, white	26091406	1
Caelo	Paraffin, white	26111203	1
Caelo	Paraffin, white	27021402	1
Bombastus	Paraffin, white	2704212549	1
Bombastus	Paraffin, white	270492	1
Bombastus	Paraffin, white	270553	2
Caelo	Paraffin, white	27061301	1
Caelo	Paraffin, white	27061401	1
Bombastus	Paraffin, white	270831	1
Caelo	Paraffin, white	27091205	1
Bombastus	Paraffin, white	270980	2
Bombastus	Paraffin, white	270981	1
Phoenix am 22.04.2014	Paraffin, white	270981	1
Bombastus	Paraffin, white	276088	2
Caelo	Paraffin, white	27121304	1
Alliance Healthcare	Paraffin, white	276088	1
Bombastus	Paraffin, white	2761127	1
Caelo	Paraffin, white	276127	1
Bombastus	Paraffin, white	276127	1
Bombastus	Paraffin, white	276968	1
Bombastus	Paraffin, white	277059	2
Gehe	Paraffin, white	277059	2
Bombastus	Paraffin, white	277060	5
Alliance Healthcare	Paraffin, white	277444	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Paraffin, white	277986	2
Phönix	Paraffin, white	277986	1
Bombastus	Paraffin, white	277987	6
Phönix	Paraffin, white	277987	1
Caelo	Paraffin, white	277987	1
Bombastus	Paraffin, white	277988	2
Gehe 02.07.2013	Paraffin, white	277988	2
Bombastus	Paraffin, white	278595	1
Gehe	Paraffin, white	279161	1
Bombastus	Paraffin, white	278620	1
Bombastus	Paraffin, white	279161	6
Bombastus	Paraffin, white	279805	1
	Paraffin, white	279969	5
Bombastus	Paraffin, white	279969	3
Gehe	Paraffin, white	279969	1
Bombastus	Paraffin, white	279970	4
AHD, 30.09.2013	Paraffin, white	279970	1
Bombastus/ Phönix	Paraffin, white	280451	1
Phönix	Paraffin, white	280451	1
Bombastus	Paraffin, white	280452	9
Caelo	Paraffin, white	280702	1
Bombastus	Paraffin, white	280702	2
Bombastus	Paraffin, white	281291	4
Phönix	Paraffin, white	281291	1
Bombastus	Paraffin, white	281539	1
Bombastus	Paraffin, white	282336	2
Bombastus	Paraffin, white	282699	1
Bombastus	Paraffin, white	282700	4
Caelo	Paraffin, white	282700	1
Bombastu-Werke AG/Gehe	Paraffin, white	282964	3
Caelo	Paraffin, white	282964	1
Bombastus	Paraffin, white	282966	7
Gehe	Paraffin, white	282966	1
Caelo	Paraffin, white	283949	1
Bombastus	Paraffin, white	283873	1
Bombastus	Paraffin, white	283948	1
Bombastus	Paraffin, white	283949	6
Bombastus	Paraffin, white	283951	1
Anzag	Paraffin, white	285231	1
Bombastus	Paraffin, white	285231	3
Bombastus	Paraffin, white	285309	1
Bombastus	Paraffin, white	285310	11
Gehe	Paraffin, white	286095	3
Bombastus	Paraffin, white	286095	4
Caelo	Paraffin, white	286095	1
Bombastus	Paraffin, white	286096	5
Gehe	Paraffin, white	286096	2
Bombastus	Paraffin, white	287149	5
Bombastus	Paraffin, white	287147	1
Euro OTC	Paraffin, white	288714	1
Bombastus	Paraffin, white	288714	4
Caelo	Paraffin, white	29051201	1
Bombastus	Paraffin, white	2909212549	1
Caelo	Paraffin, white	30012014F	1
Bombastus	Paraffin, white	30061201	1
Bombastus	Paraffin, white	3013M2549	1
Bombastus	Paraffin, white	302037	16
Bombastus	Paraffin, white	3113Q02549	1
Bombastus	Paraffin, white	4013E-02549	1
Henry Lamotte Oils GmbH	Paraffin, white	5098800	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	5098800	1
HerryLamotte/AHD	Paraffin, white	5098800	1
Bombastus	Paraffin, white	5213A2549	1
Lamotte	Paraffin, white	5215200	1
Bombastus	Paraffin, white	5215200	1
Phönix, 30.10.2012	Paraffin, white	5322600	1
Henry Lamotte Oils GmbH	Paraffin, white	5322600	1
Anzag/Lamotte	Paraffin, white	5322600	1
Lamotte	Paraffin, white	5322600	1
Henry Lamotte Oils GmbH	Paraffin, white	5489700	3
Phönix	Paraffin, white	5489700	1
gehe 26.11.13	Paraffin, white	5489700	1
Gehe am 03.02.2014	Paraffin, white	5489700	1
AHD	Paraffin, white	5489700	1
Henry Lamotte Oils GmbH	Paraffin, white	5538300	2
wollff/AHD	Paraffin, white	5538300	1
Lamotte/ Anzag	Paraffin, white	5538300	1
henry lamotte oils gmbh	Paraffin, white	5538300	1
Henry Lamotte GMBH/AHD	Paraffin, white	5691100	2
Lamotte	Paraffin, white	6000083007	1
Henry Lamotte Oils GmbH	Paraffin, white	5691100	1
Henry Lamotte Oils GmbH	Paraffin, white	6000083007	2
Bombastus	Paraffin, white	60081102	1
Caelo	Paraffin, white	90292119	1
Euro OTC	Paraffin, white	929392	1
Caelo	Paraffin, white	AKD50	1
Euro OTC	Paraffin, white	F602635	1
Euro OTC	Paraffin, white	F603290	1
Caelo	Paraffin, white	jjjjjj	1
Herbasinica	Paraffin, white		2
Phagron	Paraffin, white		1
Caelo	Paraffin, white	PV43	1
Fagron	Paraffin, white	probe	1
Caelo	Paraffin, white	4051502	1
Caelo	Paraffin, white	50115/1	1
Caelo	Paraffin, white	5121403	1
Caelo	Paraffin, white	6861393	1
Caelo	Paraffin, white	7031502	1
Caelo	Paraffin, white	11021510	1
Fagron	Paraffin, white	13/19-T05-008452	1
Caelo	Paraffin, white	133380703	1
Caelo	Paraffin, white	13341604	1
Caelo	Paraffin, white	13380704	5
cae	Paraffin, white	13380706	1
Caelo	Paraffin, white	13380706	3
Sanacorp	Paraffin, white	13380707	1
Sanacorp	Paraffin, white	13380709	1
Caelo	Paraffin, white	13380709	7
Fagron	Paraffin, white	13h19-T05-008453	1
Fagron	Paraffin, white	13119-T05-008452	1
Phönix Gotha	Paraffin, white	14087402	1
Bombastus	Paraffin, white	14087402	1
Caelo	Paraffin, white	14087406	2
Caelo	Paraffin, white	14087408	1
Caelo	Paraffin, white	14087409	3
Sanacorp	Paraffin, white	14087410	1
Euro OTC	Paraffin, white	1409009-01	3
Fagron	Paraffin, white	14124703	1
Cealo	Paraffin, white	14124707	1
Caelo	Paraffin, white	14124708	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	14189003	2
Caelo	Paraffin, white	14189007	5
Caelo	Paraffin, white	14189008	2
Fagron	Paraffin, white	14189012	2
Caelo	Paraffin, white	14189012	2
Caelo	Paraffin, white	14189019	1
Caelo	Paraffin, white	14248403	2
Caelo	Paraffin, white	14301805	1
Caelo	Paraffin, white	14324601	2
Fagron	Paraffin, white	14E30-T07	3
Fagron	Paraffin, white	14E30-T07-012094	1
Fagron	Paraffin, white	14E30-T07-012096	4
Fagron	Paraffin, white	14E30-T07-014285	8
Fagron	Paraffin, white	14j03-t003-020827	1
Euro OTC	Paraffin, white	1501037-01	4
Caelo	Paraffin, white	15050004	4
Caelo	Paraffin, white	15050005	1
Caelo	Paraffin, white	15060006	4
Caelo	Paraffin, white	15060004	6
Phönix Gotha	Paraffin, white	15060012	1
Caelo	Paraffin, white	15121407	1
Caelo	Paraffin, white	15278808	1
Caelo	Paraffin, white	15339807	11
Caelo	Paraffin, white	16181902	10
Henry Lamotte Oils	Paraffin, white	16181904	1
Caelo	Eye paraffin	17052010	1
Caelo	Paraffin, white	17101410	1
Caelo	Paraffin, white	21051509	1
Caelo	Paraffin, white	20031518	1
Caelo	Paraffin, white	22011510	1
Caelo	Paraffin, white	20061503	1
Caelo	Paraffin, white	18041504	1
Caelo	Paraffin, white	22111401	1
Caelo	Paraffin, white	25112014A	1
Caelo	Paraffin, white	28041508	1
Bombastus	Paraffin, white	287148	2
Gehe	Paraffin, white	288714	2
Bombastus	Paraffin, white	288715	3
Bombastus	Paraffin, white	288716	7
Bombastus	Paraffin, white	288717	3
Bombastus	Paraffin, white	2888717	1
Bombastus	Paraffin, white	289344	3
Phönix Gotha	Paraffin, white	289344	1
Bombastus	Paraffin, white	290271	4
Fagron	Paraffin, white	290441	1
Bombastus	Paraffin, white	290441	4
Caelo	Paraffin, white	290441	2
Caesar & Loretz GmbH	Paraffin, white	290441	2
Bombastus	Paraffin, white	290714	3
Bombastus	Paraffin, white	291318	14
Bombastus	Paraffin, white	292219	2
Bombastus	Paraffin, white	292220	3
Bombastus	Paraffin, white	295301	6
Bombastus	Paraffin, white	297282	5
Bombastus	Paraffin, white	300120	9
Caelo	Paraffin, white	30121427	1
Euro OTC	Paraffin, white	5691100	1
Sanacorp	Paraffin, white		1
Caelo	Paraffin, white	323915339805	1
Fagron	Paraffin, white	27930	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	4091517	1
Caelo	Paraffin, yellow	6895914	1
Sanacorp	Paraffin, white	106/AKD4	1
Bombastus	Paraffin, yellow	106429	1
Sanacorp	Paraffin, white	107/AKD6	1
Caelo	Paraffin, white	108/AKD12	1
Fagron	Paraffin, white	1166979	1
Caelo	Paraffin, white	1111111	1
Bombastus	Paraffin, white	111967	1
Bombastus	Paraffin, white	12/2019	1
Caelo	Paraffin, white	12081504	1
Caelo	Paraffin, yellow	12306302	1
Caelo	Paraffin, yellow	12321902	1
Euro OTC	Paraffin, yellow	12321908	1
Caelo	Paraffin, yellow	13365104	1
Gehe 20.06.15/euro	Paraffin, white	13380709	1
Noweda	Paraffin, white	13380721	1
Fagron	Paraffin, yellow	13K20-T05-009125	1
Caelo	Paraffin, white	14092015B	1
Caelo	Paraffin, yellow	14248403	3
Caelo	Paraffin, yellow	14324301	6
Caesar&Loretz/Noweda	Paraffin, yellow	14324301	1
Phönix	Paraffin, yellow	14324301	1
Caelo	Paraffin, yellow	14324402	1
Caelo	Paraffin, yellow	14324501	10
Caesar & Loretz GmbH Hi...	Paraffin, yellow	14324601	1
Caelo	Paraffin, yellow	14324601	3
Caelo	Paraffin, yellow	1432501	1
Fagron	Paraffin, white	14F30-T07-014285	1
Fagron	Paraffin, white	14I03-T03-020831	1
Fagron	Paraffin, white	14J03-020827	1
Fagron	Paraffin, white	14j03-t03-020827	1
Fagron	Paraffin, white	14J03-T03	1
Fagron	Paraffin, white	14j03-t03020831	1
Fagron	Paraffin, white	14J03T-03-020831	1
Fagron	Paraffin, white	14J03-T03-020831	2
Fagron	Paraffin, white	14J03T03020831	1
Fragon/Phoenix	Paraffin, white	14J03-T03-020831	1
Fagron	Paraffin, yellow	14K27-T29-019303	2
Fagron	Paraffin, white	14L11-T03-027930	6
Eurho/Noweda	Paraffin, white	1501037-01	1
Euro OTC	Paraffin, yellow	1503036-01	6
Caelo	Paraffin, yellow	1503036-01	1
Euro OTC	Paraffin, white	1505017-01	5
Hedinger	Paraffin, white	15060004	1
Caelo	Paraffin, white	15060007	2
Caelo	Paraffin, white	15060011	2
Gehe	Paraffin, white	15060012	1
Anzag/Henry Lamotte	Paraffin, white	15060014	1
Caelo	Paraffin, yellow	15063902	2
Caelo	Paraffin, yellow	15063904	1
Sana 29.11.Ek.6,95	Paraffin, yellow	15063904	1
Caelo	Paraffin, yellow	150717VF	1
Euro OTC	Paraffin, white	1509019-01	6
Euro OTC	Paraffin, white	150901901	1
Ebert & Jacobi 19.07...	Paraffin, white	1509019-01	1
Caelo	Paraffin, white	15162501	2
Caelo	Paraffin, white	15162506	1
Caelo	Paraffin, white	15162519	1
Sanacorp	Paraffin, white	15204301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, yellow	15219903	5
Caelo	Paraffin, yellow	15219904	1
Caelo	Paraffin, white	15278809	2
Caelo	Paraffin, yellow	15330401	8
Caesar & Loretz/Noweda	Paraffin, yellow	15330401	1
Bombastus	Paraffin, yellow	15330401	1
Caelo	Paraffin, yellow	15330402	4
Noweda	Paraffin, yellow	15330402	1
Caelo	Paraffin, yellow	15330502	1
Caelo	Paraffin, yellow	15330601	12
WE:Sanacorp 05.09.2016...	Paraffin, yellow	15330601	2
Caelo	Paraffin, white	15336808	2
Sana 21.03.16 EK 4,96	Paraffin, white	15336810	1
Caelo	Paraffin, white	15339801	1
Caelo	Paraffin, white	15339805	3
Noweda	Paraffin, white	15339806	1
Euro OTC	Paraffin, white	15339806	1
gehe	Paraffin, white	15339806	1
Caesar & Loretz/Noweda	Paraffin, white	15339806	1
Caelo	Paraffin, white	15339808	12
Caesar & Loretz	Paraffin, white	15339808	1
Caelo	Paraffin, white	15349702	7
Caelo	Paraffin, white	15349705	1
Caelo	Paraffin, white	15349706	4
Phoenix	Paraffin, white	15349706	1
Caelo	Paraffin, white	15349707	4
Caelo	Paraffin, white	15349709	3
Caelo	Paraffin, white	15359105	1
Caelo	Paraffin, white	15355806	1
Sanacorp	Paraffin, white	15441504	1
Caelo	Paraffin, white	15441504	2
Caelo	Paraffin, white	15441505	2
Caelo	Paraffin, white	15441509	1
Caelo	Paraffin, white	15441506	1
Caelo	Paraffin, white	15441511	3
Caelo	Paraffin, white	15441515	1
Caelo	Paraffin, white	1544159	1
Euro OTC	Paraffin, white	1602009-01	3
Euroo/Gehe	Paraffin, white	1602009-01	1
Caelo	Paraffin, white	16029210	2
Phoenix	Paraffin, white	16029210	1
Caelo	Paraffin, white	16029208	1
Caelo	Paraffin, white	16029212	1
Caelo	Paraffin, white	16029219	4
Caelo	Paraffin, white	16029220	4
Caelo	Paraffin, white	16029222	3
Caelo	Paraffin, yellow	16032502	2
Noweda	Paraffin, yellow	16032502	1
Phönix	Paraffin, yellow	16032602	1
Caelo	Paraffin, yellow	16032602	5
Caelo	Paraffin, yellow	16032603	3
Sanacorp, 21.07.2016	Paraffin, white	1605019-01	2
Euro OTC	Paraffin, white	1605019-01	3
Euro OTC	Paraffin, white	160519-01	1
Caesar & Loretz/ Noweda	Paraffin, yellow	1606913	1
Euro OTC	Paraffin, white	1608041-01	3
Caelo	Paraffin, white	16082701	2
Caelo	Paraffin, yellow	16096903	1
Ichthyol-Gesellschaft	Paraffin, white	160829VA	1
Caelo	Paraffin, yellow	16096904	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar&Loretz	Paraffin, yellow	16096903	1
Caelo	Paraffin, yellow	16096913	1
Euro OTC	Paraffin, white	1612007-01	1
Caelo	Paraffin, white	161338003	1
Ichthyol/Phönix	Paraffin, white	16138003	1
phönix	Paraffin, white	16138003	1
Noweda	Paraffin, white	16138004	1
Caelo	Paraffin, white	16138004	13
Caelo	Paraffin, yellow	161587801	1
Caelo	Paraffin, white	16161901	1
Caelo	Paraffin, white	16181901	5
Celo	Paraffin, white	16181902	2
Alte Apotheke	Paraffin, white	16181902	1
Caelo	Paraffin, white	16181903	6
Caelo	Paraffin, white	16181904	10
Allianz Healthcare	Paraffin, white	16181904	1
Caelo	Paraffin, yellow	1618701	1
Caelo	Paraffin, yellow	16187702	4
Caelo	Paraffin, yellow	16187801	6
Caelo	Paraffin, white	16227401	2
Caelo	Paraffin, white	16227402	6
PHÖNIX	Paraffin, white	16227402	1
Caelo	Paraffin, white	16227403	4
Caesar & Loretz GmbH Hi...	Paraffin, white	16227403	1
Caelo	Paraffin, white	16227404	12
Caelo	Paraffin, white	16227405	3
Caesar&Loretz GmbH,27.0...	Paraffin, white	16227405	1
Caelo	Paraffin, white	16227406	7
Caelo	Paraffin, white	16227411	4
Fagron	Paraffin, white	162274112	4
Caelo	Paraffin, white	16227412	9
Sanacorp	Paraffin, white	16227412	1
Caelo	Paraffin, white	16241105	2
Caelo	Paraffin, white	16271807	1
Caelo	Paraffin, white	16271810	2
Caelo	Paraffin, white	16272001	1
Caelo	Paraffin, white	16272002	2
Sanacorp 37.03.2017 ek:...	Paraffin, white	16272002	1
Caelo	Paraffin, white	16272003	2
Euro OTC	Paraffin, white	16272003	1
Caelo	Paraffin, white	16272004	3
Fagron	Paraffin, white	16272004	3
Cealo/AHD	Paraffin, white	16272004	1
Caelo	Paraffin, white	16272008	2
Caelo	Paraffin, white	16272017	11
Caelo	Paraffin, white	16272010	1
Caelo	Paraffin, yellow	16281001	10
Gehe	Paraffin, yellow	16281001	1
Caelo	Paraffin, yellow	16281002	1
Caelo	Paraffin, white	16326702	10
Caesar & Loretz/Noweda	Paraffin, white	16326702	1
Caelo	Paraffin, white	16326707	1
Fagron	Paraffin, yellow	16E17-T01-035266	1
Fagron	Paraffin, white	16L27-T01-042281	2
Caelo	Paraffin, white	1672017	1
Caelo	Paraffin, white	170203VA/16181901	1
Caelo	Paraffin, white	17038201	1
Caelo	Paraffin, white	17038202	3
Caelo	Paraffin, white	17038203	2
Caelo	Paraffin, white	17038206	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, white	17038204	1
Caelo	Paraffin, white	17038213	1
Caelo	Paraffin, white	17038215	2
Fagron	Paraffin, white	17052001	1
Caelo	Paraffin, white	17052010	1
Caelo	Paraffin, yellow	17052401	3
Caelo	Paraffin, yellow	17052402	1
Caelo	Paraffin, white	17106001	1
Caesar & Loretz/ Noweda	Paraffin, white	1710005	1
Caelo	Paraffin, yellow	17082015E	1
chaelo/gehe	Paraffin, white	17106002	1
Caelo	Paraffin, white	17106002	1
Caelo	Paraffin, white	171114VA/16241105	1
Caelo	Paraffin, yellow	17130402	2
Caelo	Paraffin, white	1720401	1
Caelo	Paraffin, white	17204201	1
Caelo	Paraffin, white	17204203	1
Caelo	Paraffin, yellow	17223102	3
GEHE	Paraffin, yellow	17223102	1
Caelo	Paraffin, white	17347901	2
Caelo	Paraffin, white	17350003	1
Caelo	Paraffin, white	17350004	1
Caelo	Paraffin, white	20071505	1
Caelo	Paraffin, yellow	28071509	1
Bombastus	Paraffin, yellow	28270	1
Bombastus	Paraffin, yellow	285309	1
Bombastus	Paraffin, white	282964	2
Bombastus	Paraffin, yellow	285310	1
Euro OTC	Paraffin, yellow	289344	1
Bombastus	Paraffin, yellow	289344	1
Bombastus	Paraffin, yellow	290713	1
Bombastus	Paraffin, yellow	290714	1
Fagron	Paraffin, white	291318	2
Bombastus	Paraffin, white	291319	1
Caelo	Paraffin, white	292220	1
Bombastus	Paraffin, white	293553	1
Bombastus	Paraffin, white		1
Bombastus	Paraffin, yellow	293629	11
Bombastus	Paraffin, white	294852	4
Bombastus	Paraffin, white	294851	1
Caelo	Paraffin, white	2948852	1
Bombastus	Paraffin, white	295265	2
Bambustus/Noweda	Paraffin, white	295265	1
Caelo	Paraffin, white	295300	1
Bombastus	Paraffin, white	295344	4
Bombastus	Paraffin, yellow	296654	9
Bambustus/Noweda	Paraffin, white	296854	4
Caelo	Paraffin, white	297283	1
Bombastus	Paraffin, white	297283	7
Euro OTC	Paraffin, white	297283	1
Bombastus	Paraffin, white	298427	2
Bombastus	Paraffin, white	298599	10
Bombastus	Paraffin, white	298600	1
Bombastus	Paraffin, white	299468	2
Bombastus	Paraffin, white	299469	5
Caelo	Paraffin, white	299469	1
Bombastus	Paraffin, white	2999469	1
Bombastus	Paraffin, white	300119	3
GEHE	Paraffin, white	300119	1
Bombastus	Paraffin, yellow	300422	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Paraffin, yellow	300423	6
Caelo	Paraffin, white	301424	1
Bombastus	Paraffin, yellow	301908	6
WE:Phönix 02.09.2016 E...	Paraffin, white	302037	2
Bambynek/Gehe	Paraffin, white	303786	1
Caelo	Paraffin, white	303786	1
Bombastus	Paraffin, yellow	303809	3
Bombastus	Paraffin, white	304124	8
Bombastus	Paraffin, white	304125	1
Bombastus	Paraffin, yellow	304847	4
Bombastus	Paraffin, white	305140	1
Fagron	Paraffin, white	306019	2
Bombastus	Paraffin, white	306020	6
Caelo	Paraffin, white	306020	1
apomix/Noweda	Paraffin, white	306020	1
Bombastus	Paraffin, yellow	306182	2
Bombastus	Paraffin, yellow	307079	5
Bombastus	Paraffin, white	307661	6
Caelo	Paraffin, white	307661	1
Bomastus/Phönix	Paraffin, white	307661	1
Noweda	Paraffin, white	307661	1
Bombastus	Paraffin, white	308197	2
Bombastus	Paraffin, white	309563	2
Bombastus	Paraffin, white	309564	2
Sanacorp E.K:5,64 -100...	Paraffin, white	309564	1
Gehe	Paraffin, yellow	310191	1
Lamotte 001	Paraffin, white	5691100	1
Sanacorp	Paraffin, yellow	57/AKD4	1
Phoenix	Paraffin, yellow	57/AKD7	1
Sanacorp	Paraffin, yellow	58/AKD12	1
Caelo	Paraffin, white	6000083007	2
Lamotte	Paraffin, white	6002062001	2
HenryLamotte	Paraffin, white	6002062001	3
Caelo	Paraffin, white	6002062001	1
Gehe	Paraffin, white	60081104	1
ACM	Paraffin, white	KNC1628J15	1
Caelo	Paraffin, white	14201003	1
Beiersdorf AG	Eye paraffin	15441502	1
sanacorp, 6,08EUR, 29.08...	Paraffin, white	16138003	1
Sanacorp,04.04.17, 6,08EUR	Paraffin, white	16181904	1
Bombastus	Paraffin, white	267671	2
GEHE	Paraffin, white	270163	1
GEHE	Paraffin, white	270834	1
Bombastus	Paraffin, white	270984	1
Caelo	Paraffin, white	12081502	1
Caelo	Paraffin, white	1210006-01	1
Gehe	Paraffin, white	1210006-01	1
Caelo	Paraffin, white	1215106	2
Sanacorp	Paraffin, white	12374502	1
Gutenberg-Apotheke	Paraffin, white	13185502	1
VDL;13.06.14	Paraffin, white	13268421	1
Gehe	Paraffin, white	277060	1
Caelo	Paraffin, white	51210195	1
Henry Lamotte/N/14.12.12	Paraffin, white	5322600	1
Fagron	Paraffin, white	790619-013771	1
Caelo	Paraffin, white	P159/14	1
GEHE	Paraffin, liquid	12102101	1
Ichtyol, Noweda	Unguentum Cordes®	12A147	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	414610	1
Bombastus	Wool alcohol ointment DAB	274282	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sonnen-Apotheke Wefensl...	Wool alcohol ointment DAB	1/16	1
Apomix	Wool alcohol ointment DAB	1A170214	1
Apomix/Noweda	Wool alcohol ointment SR wit...	1A180215	1
apomix	Wool alcohol ointment DAB	1A180215	1
Apomix	Wool alcohol ointment SR wit...	2A070515	1
Caelo	Wool alcohol ointment DAB	4051505	1
Caelo	Wool alcohol ointment DAB	4067201	1
Apomix	Wool alcohol ointment SR wit...	5231213	1
Apomix	Wool alcohol ointment SR wit...	6A290714	1
Henry Lamotte Oels GmbH	Wool alcohol ointment SR wit...	7A111214	2
Sonnen-Apotheke Wefensl...	Wool alcohol ointment DAB	9/17	1
Defektur	Wool alcohol ointment DAB	10/15	1
Bombastus	Wool alcohol ointment DAB	104374	1
Bombastus	Wool alcohol ointment SR wit...	106236	1
Bombastus	Wool alcohol ointment SR/DAC	108187	1
Beiersdorf AG	Wool alcohol ointment DAB	11062014A	1
Euro OTC	Wool alcohol ointment DAB	1204030-01	1
Eurho/Noweda	Wool alcohol ointment DAB	1211017-01	1
Bombastus	Wool alcohol ointment DAB	1211017-01	1
Caelo	Wool alcohol ointment DAB	13067103	7
Ebert	Wool alcohol ointment DAB	13242602	1
Caelo	Wool alcohol ointment DAB	13242902	6
Caelo	Wool alcohol ointment DAB	1324602	1
Anzag 25.02.2014 Ek 8...	Wool alcohol ointment DAB	13255001	1
Caelo	Wool alcohol ointment DAB	1325501	1
DAB 2008	Wool alcohol ointment DAB	13255101	1
Caelo	Wool alcohol ointment SR wit...	13264102	2
Caele/ Jenne	Wool alcohol ointment SR wit...	13340003	1
Fagron	Wool alcohol ointment DAB	13D09-T04	1
Caelo	Wool alcohol ointment SR wit...	14032502	1
Caelo	Wool alcohol ointment DAB	14022902	1
Caelo	Wool alcohol ointment DAB	14067001	5
Phönix	Wool alcohol ointment DAB	14067001	1
Caelo	Wool alcohol ointment DAB	1406701	1
Caelo	Wool alcohol ointment DAB	14067101	5
Berkel AHK	Wool alcohol ointment DAB	14067201	1
Caelo	Wool alcohol ointment DAB	14067301	1
Caelo	Wool alcohol ointment DAB	14067401	4
Caelo	Wool alcohol ointment DAB	14152901	1
Caelo	Wool alcohol ointment DAB	14152803	1
Caelo	Wool alcohol ointment DAB	14152902	1
Phoenix	Wool alcohol ointment DAB	14153101	1
Sana 13.11.14 EK 8,96	Wool alcohol ointment DAB	14153101	1
Caelo	Wool alcohol ointment DAB	14153601	8
Fagron	Wool alcohol ointment DAB	14153601	1
Caelo	Wool alcohol ointment DAB	14207702	1
Caelo	Wool alcohol ointment DAB	14283001	14
Caesar & Loretz GmbH Hi...	Wool alcohol ointment DAB	14283001	1
Caelo	Wool alcohol ointment DAB	14283102	9
Caelo	Wool alcohol ointment DAB	14283103	11
Caelo	Wool alcohol ointment DAB	1428801	1
Caelo	Wool alcohol ointment DAB	143020402	1
Fagron	Wool alcohol ointment DAB	14H19-T01-011354	1
Alliance Health 30.10.2...	Wool alcohol ointment DAB	14H19-T02-011354	1
Fagron	Wool alcohol ointment DAB	14H19-T02-011354	1
Fagron	Wool alcohol ointment DAB	14H19-T02-018184	2
Fagron	Wool alcohol ointment DAB	14H19T02018184	1
Caelo	Wool alcohol ointment DAB	15035501	5
Caelo	Wool alcohol ointment SR wit...	15047002	2
Caelo	Wool alcohol ointment SR wit...	15047004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol ointment SR wit...	15047003	1
Caelo	Wool alcohol ointment SR wit...	15047101	6
Caelo	Wool alcohol ointment DAB	15225501	6
P.W. Beyvers GmbH/Phönix	Wool alcohol ointment DAB	15225501	3
Phönix Gotha	Wool alcohol ointment DAB	15225601	1
caesar&lorenz,Phönix	Wool alcohol ointment DAB	15225601	1
Fagron	Wool alcohol ointment DAB	15225601	1
Caelo	Wool alcohol ointment DAB	1522802	1
Sanacorp 31.10.15	Wool alcohol ointment DAB	15228202	2
Caelo	Wool alcohol ointment SR wit...	15334703	1
Noweda	Wool alcohol ointment DAB	15361801	1
Caelo	Wool alcohol ointment DAB	15361901	8
Caeolo/Noweda	Wool alcohol ointment DAB	15361901	1
Caelo	Wool alcohol ointment DAB	15362001	10
Caelo	Wool alcohol ointment DAB	15362002	8
Caesar&Loretz GmbH	Wool alcohol ointment DAB	15362002	1
Noweda	Wool alcohol ointment DAB	15362002	1
Caelo	Wool alcohol ointment DAB	15362003	7
Caelo	Wool alcohol ointment DAB	154035502	1
Caelo	Wool alcohol ointment DAB	15407001	1
Fagron	Wool alcohol ointment DAB	15E18-20-021463	1
Fagron	Wool alcohol ointment DAB	15E18-T20-021463	2
Caelo	Wool alcohol ointment DAB	16007202	6
Caelo	Wool alcohol ointment DAB	16007203	1
Caelo	Wool alcohol ointment DAB	16066501	7
Caelo	Wool alcohol ointment DAB	16066401	1
Caelo	Wool alcohol ointment DAB	16066701	4
Caelo	Wool alcohol ointment DAB	16066801	7
Euro OTC	Wool alcohol ointment DAB	1608040-01	1
Caelo	Wool alcohol ointment SR wit...	16142802	2
Caelo	Wool alcohol ointment SR wit...	16142803	1
Caelo	Wool alcohol ointment DAB	16177603	1
Caelo	Wool alcohol ointment DAB	16198101	9
Caesar & Loretz	Wool alcohol ointment DAB	16198401	1
Caelo	Wool alcohol ointment DAB	16198401	4
AHD 06.02.17 EK.9,28	Wool alcohol ointment DAB	16198401	1
Caelo	Wool alcohol ointment SR/DAC	16198401	1
Sanacorp	Wool alcohol ointment DAB	16198401	1
Caelo	Wool alcohol ointment DAB	16198501	12
Caealo	Wool alcohol ointment DAB	16198501	1
phönix	Wool alcohol ointment DAB	16198501	1
Bombastus	Wool alcohol ointment DAB	16198502	1
Caelo	Wool alcohol ointment DAB	16198601	7
Gehe	Wool alcohol ointment DAB	16198601	1
Caelo	Wool alcohol ointment DAB	16248802	2
Caelo	Wool alcohol ointment SR wit...	1625102	1
Caelo	Wool alcohol ointment SR wit...	16251302	1
Caelo	Wool alcohol ointment DAB	16295101	6
Caelo	Wool alcohol ointment DAB	16295102	3
Caesar & Loretz/Noweda	Wool alcohol ointment DAB	16295102	1
Caelo	Wool alcohol ointment DAB	16298502	1
GEHE	Wool alcohol ointment DAB	1633170	1
phönix	Wool alcohol ointment DAB	16331701	2
Caelo	Wool alcohol ointment DAB	16331701	2
Caelo/Gehe	Wool alcohol ointment DAB	16331704	2
Caelo	Wool alcohol ointment DAB	16331704	1
Caelo	Wool alcohol ointment DAB	16331802	1
Caelo	Wool alcohol ointment DAB	16331901	5
Fargon	Wool alcohol ointment DAB	16H03-T01-035270	1
Bombastus	Wool alcohol ointment DAB	17001875	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol ointment DAB	1704903	1
Caelo	Wool alcohol ointment DAB	17104901	3
Phönix	Wool alcohol ointment DAB	17104901	2
Sanacoro W.E:22.02.18 ...	Wool alcohol ointment DAB	17104902	1
Caelo	Wool alcohol ointment DAB	17104903	1
Caleo	Wool alcohol ointment DAB	17104903	1
Caelo	Wool alcohol ointment DAB	17105001	4
Caesar & Loretz GmbH, P...	Wool alcohol ointment DAB	17105001	1
Caelo	Wool alcohol ointment DAB	17105401	3
GEHE	Wool alcohol ointment DAB	17141202	1
Caelo	Wool alcohol ointment DAB	17141203	4
Caelo	Wool alcohol ointment DAB	17141202	1
Caelo	Wool alcohol ointment DAB	17141401	2
Caelo	Wool alcohol ointment DAB	17141404	1
Caelo	Wool alcohol ointment DAB	17141601	1
Defektur	Wool alcohol ointment DAB	18/14	1
Bombastus	Wool alcohol ointment DAB	18001096	1
Caelo	Wool alcohol ointment DAB	21051510	1
GEHE 17.06.2014	Wool alcohol ointment DAB	23033310	1
Pharmachem	Wool alcohol ointment SR wit...	250116	1
Caelo	Wool alcohol ointment DAB	259166	1
Defektur	Wool alcohol ointment DAB	26/14	1
Bombastus	Wool alcohol ointment DAB	277024	2
Bombastus	Wool alcohol ointment SR wit...	277024	4
Gehe	Wool alcohol ointment SR wit...	279105	2
Bombastus	Wool alcohol ointment SR wit...	280057	1
Bombastus	Wool alcohol ointment DAB	280061	1
Bombastus	Wool alcohol ointment SR wit...	283174	1
Bombastus	Wool alcohol ointment DAB	283176	1
Bombastus	Wool alcohol ointment DAB	283618	8
gehe	Wool alcohol ointment DAB	284378	1
Bombastus	Wool alcohol ointment SR wit...	284378	2
Bombastus	Wool alcohol ointment SR wit...	284379	1
Phönix Gotha	Wool alcohol ointment SR wit...	284381	1
Bombastus	Wool alcohol ointment SR wit...	284381	1
Bombastus	Wool alcohol ointment DAB	285423	1
Apomix	Wool alcohol ointment DAB	285950	1
Bombastus	Wool alcohol ointment SR wit...	286880	1
Caelo	Wool alcohol ointment DAB	286882	1
Bombastus	Wool alcohol ointment SR wit...	286884	1
Phönix Gotha	Wool alcohol ointment DAB	286884	1
Bombastus	Wool alcohol ointment DAB	287013	4
Noweda	Wool alcohol ointment DAB	287013	1
Bombastus	Wool alcohol ointment SR wit...	288686	1
Gehe	Wool alcohol ointment SR wit...	288687	1
Bombastus	Wool alcohol ointment SR wit...	288688	1
Phoenix 27.11.14 EK 6,08	Wool alcohol ointment SR wit...	288688	1
Phönix Gotha	Wool alcohol ointment DAB	288688	1
Bombastus	Wool alcohol ointment DAB	288688	1
Sanacop	Wool alcohol ointment SR wit...	288688	1
Bomba/ Kehr	Wool alcohol ointment SR wit...	288689	1
Bombastus	Wool alcohol ointment DAB	288766	1
Bombastus	Wool alcohol ointment DAB	288834	8
UNISAN	Wool alcohol ointment DAB	289509	1
Bombastus	Wool alcohol ointment DAB	289509	6
Sanacorp	Wool alcohol ointment DAB	289509	1
Gehe	Wool alcohol ointment DAB	289510	1
	Wool alcohol ointment SR wit...	291175	2
Bombastus	Wool alcohol ointment SR wit...	291175	1
Caelo	Wool alcohol ointment DAB	291175	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol ointment SR wit...	291176	1
Bombastus	Wool alcohol ointment DAB	291794	1
Bombastus	Wool alcohol ointment DAB	291796	3
Bombastus	Wool alcohol ointment DAB	291797	2
Bombastus	Wool alcohol ointment DAB	293091	2
Bombastus	Wool alcohol ointment DAB	293092	6
Caesar und Loretz	Wool alcohol ointment DAB	293092	2
Caelo	Wool alcohol ointment DAB	293093	1
Bombastus	Wool alcohol ointment DAB	29309932ma1	1
Bombastus	Wool alcohol ointment SR wit...	293374	1
Bombastus	Wool alcohol ointment DAB	293376	2
Bombastus	Wool alcohol ointment SR wit...	293378	1
Bombastus	Wool alcohol ointment DAB	293956	10
Bombastus	Wool alcohol ointment DAB	295357	1
Bombastus	Wool alcohol ointment SR wit...	295358	2
Bombastus	Wool alcohol ointment DAB	295358	1
Bombastus	Wool alcohol ointment SR wit...	295359	1
Bombastus	Wool alcohol ointment SR/DAC	295359	1
Bombastus	Wool alcohol ointment SR wit...	295360	2
Bombastus	Wool alcohol ointment DAB	296235	1
Bombastus	Wool alcohol ointment DAB	296435	8
Bombastus	Wool alcohol ointment DAB	296436	3
Bombastus	Wool alcohol ointment SR wit...	297125	1
Bombastus	Wool alcohol ointment SR wit...	297127	1
Bombastus	Wool alcohol ointment SR wit...	297128	2
Bombastus	Wool alcohol ointment SR/DAC	297128	1
Bombastus	Wool alcohol ointment DAB	297129	1
Bombastus	Wool alcohol ointment SR wit...	297130	1
Bombastus	Wool alcohol ointment DAB	298577	9
Bombastus	Wool alcohol ointment DAB	298578	1
Bombastus	Wool alcohol ointment SR/DAC	300080	1
Bombastus	Wool alcohol ointment SR wit...	300083	2
Bombastus	Wool alcohol ointment SR/DAC	300082	1
Bombastus	Wool alcohol ointment SR/DAC	300084	3
Bombastus	Wool alcohol ointment DAB	300281	1
Bombastus	Wool alcohol ointment DAB	300282	5
Caelo	Wool alcohol ointment DAB	300283	1
Bombastus	Wool alcohol ointment DAB	300283	3
Bombastus	Wool alcohol ointment SR wit...	30079	1
Bombastus	Wool alcohol ointment DAB	301200	2
Bombastus	Wool alcohol ointment DAB	301199	1
Bombastus	Wool alcohol ointment SR/DAC	301935	1
Bombastus	Wool alcohol ointment SR/DAC	301937	1
Bombastus	Wool alcohol ointment SR/DAC	301940	1
Bombastus	Wool alcohol ointment SR/DAC	301938	1
Bombastus	Wool alcohol ointment DAB	302206	9
Bombastus	Wool alcohol ointment DAB	302208	1
Bombastus	Wool alcohol ointment DAB	30226	2
Bombastus	Wool alcohol ointment SR/DAC	302598	2
Apomix/Gehe	Wool alcohol ointment SR wit...	302599	3
Bombastus	Wool alcohol ointment SR/DAC	302599	3
Bombastus	Wool alcohol ointment SR/DAC	302600	1
Bombastus	Wool alcohol ointment DAB	303846	11
Bombastus	Wool alcohol ointment DAB	303847	2
Bombastus	Wool alcohol ointment DAB	303848	2
Bombastus	Wool alcohol ointment SR/DAC	303993	1
Bombastus	Wool alcohol ointment SR/DAC	303994	1
Bombastus	Wool alcohol ointment SR/DAC	303995	3
Bombastus	Wool alcohol ointment SR/DAC	303996	1
Bombastus	Wool alcohol ointment DAB	304998	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol ointment SR/DAC	305546	1
Bombastus	Wool alcohol ointment SR/DAC	305548	1
Bombastus	Wool alcohol ointment DAB	305878	1
Bombastus	Wool alcohol ointment DAB	305879	13
Bombastus	Wool alcohol ointment DAB	305880	3
Caelo	Wool alcohol ointment DAB	305881	1
Bombastus	Wool alcohol ointment DAB	307774	2
Bombastus	Wool alcohol ointment DAB	307775	2
Bombastus	Wool alcohol ointment DAB	307776	1
Caelo	Wool alcohol ointment SR/DAC	308141	1
Bombastus	Wool alcohol ointment SR/DAC	308144	2
Bombastus	Wool alcohol ointment SR/DAC	308145	1
Bombastus	Wool alcohol ointment SR wit. ...	308145	1
Bombastus	Wool alcohol ointment SR/DAC	308146	1
GEHE	Wool alcohol ointment SR/DAC	310488	2
Bombastus	Wool alcohol ointment DAB	311388	3
Beiersdorf/ Sanacorp	Wool alcohol ointment DAB	34344210	1
Noweda / Beiersdorf AG	Wool alcohol ointment DAB	34424610	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	34554710	1
Noweda	Wool alcohol ointment DAB	34554710	1
Defektur	Wool alcohol ointment DAB	38/14	1
Beiersdorf Gehe	Wool alcohol ointment DAB	3841	1
Apotheke/Defekturabfüll. ...	Wool alcohol ointment DAB		1
Beiersdorf AG/Geilenkir. ...	Wool alcohol ointment DAB	40726610	1
Baiersdorf AG	Wool alcohol ointment DAB	40726610	1
Gehe/Beiersdorf	Wool alcohol ointment DAB	40726610	1
Beiersdorf/ Gehe	Wool alcohol ointment DAB	40916810	1
Fagron	Wool alcohol ointment DAB	40956810	1
Defektur	Wool alcohol ointment DAB	41/14	1
Caelo	Wool alcohol ointment DAB	41910710	1
Caelo	Wool alcohol ointment DAB	41047010	1
Beiersdorf	Wool alcohol ointment DAB	41910710	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	41910710	1
Caesar & Loretz GmbH	Wool alcohol ointment DAB	41930610	1
Gehe/Beiersdorf	Wool alcohol ointment DAB	41930610	2
Phönix Gotha	Wool alcohol ointment DAB	41930610	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	41940610	2
Caelo	Wool alcohol ointment DAB	41940610	1
Phoenix	Wool alcohol ointment DAB	41940610	1
Beiersdorf / E	Wool alcohol ointment DAB	41940610	1
Gehe	Wool alcohol ointment DAB	41950710	1
Defektur	Wool alcohol ointment DAB	42/14	1
Beiersdorf/Anzag	Wool alcohol ointment DAB	42531910	1
Sanacorp	Wool alcohol ointment DAB	44/AKD10	1
Beiersdorf/ Sanacorp	Wool alcohol ointment DAB	44824310	1
Defektur	Wool alcohol ointment DAB	45/14	1
Beiersdorf/ Nowea	Wool alcohol ointment DAB	45054410	1
Defektur	Wool alcohol ointment DAB	5/15	1
Caelo	Wool alcohol ointment DAB	50/AKD7	1
Gehe	Wool alcohol ointment DAB	50444710	1
Beiersdorf/ Gehe	Wool alcohol ointment DAB	50414810	1
Gehe/Beiersdorf	Wool alcohol ointment DAB	50444710	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	50444710	1
Beiersdorf	Wool alcohol ointment DAB	50454710	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	50815310	1
Beiersdorf/Beiersdorf	Wool alcohol ointment DAB	50815310	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	50815310	1
Gehe/Beiersdorf	Wool alcohol ointment DAB	50815310	1
Beiersdorf / Fiebig	Wool alcohol ointment DAB	50815310	1
Caelo	Wool alcohol ointment DAB	50835210	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Beiersdorf AG	Wool alcohol ointment DAB	50845210	1
Beiersdorf/N	Wool alcohol ointment DAB	50845210	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	50845210	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	50845210	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	50935410	2
Beiersdorf/Phoenix	Wool alcohol ointment DAB	50935410	1
Beiersdorf/N	Wool alcohol ointment DAB	51015610	1
Beiersdorf über Phoenix	Wool alcohol ointment DAB	53437010	1
Fagron	Wool alcohol ointment DAB	53437010	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	54127710	1
Bayersdorf	Wool alcohol ointment DAB	54137710	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	54157810	1
Henry Lamotte Oils GmbH	Wool alcohol ointment DAB	5860600	1
Bombastus	Wool alcohol ointment DAB	60080202	2
Bombastus	Wool alcohol ointment DAB	60080203	2
Bombastus	Wool alcohol ointment SR/DAC	60089603	1
Beiersdorf AG	Wool alcohol ointment DAB	60419110	1
Caelo	Wool alcohol ointment DAB	60429110	1
Gehe/Beiersdorf	Wool alcohol ointment DAB	60429110	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	60539210	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	60539210	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	61149910	1
Phönix	Wool alcohol ointment DAB	61149910	1
Beiersdorf/AHD	Wool alcohol ointment DAB	61149910	1
Beiersdorf/ Noweda	Wool alcohol ointment DAB	61841610	1
Beiersdorf/Noweda	Wool alcohol ointment DAB	62922210	1
PHÖNIX	Wool alcohol ointment DAB	63052310	1
Beiersdorf	Wool alcohol ointment DAB	64332810	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	64342810	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	64342810	1
Phönix	Wool alcohol ointment DAB	64833610	1
Medphano Arzneimittel /...	Wool alcohol ointment SR wit...	6550115	1
medphano Arzneimittel	Wool alcohol ointment SR/DAC	6550115	1
Beiersdorf	Wool alcohol ointment DAB	70314010	1
Caelo	Wool alcohol ointment DAB	70314010	2
Euro OTC	Wool alcohol ointment DAB	70434110	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	71114510	1
beiersdorf/ Noweda	Wool alcohol ointment DAB	71214810	1
Beiersdorf/GEHE	Wool alcohol ointment DAB	72915810	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	74017010	1
Beiersdorf	Wool alcohol ointment DAB	74747510	1
Beiersdorf/Sanacorp	Wool alcohol ointment DAB	74757510	1
Beiersdorf	Wool alcohol ointment DAB	74757510	1
Beiersdorf / AEP	Wool alcohol ointment DAB	81819410	1
Beiersdorf/ Phoenix	Wool alcohol ointment DAB	81819410	1
Bombastus	Wool alcohol ointment DAB	251850	1
Bombastus	Wool alcohol ointment DAB	254685	1
Bombastus	Wool alcohol ointment DAB	252186	1
Bombastus	Wool alcohol ointment DAB	256518	1
Bombastus	Wool alcohol ointment DAB	261436	1
Bombastus	Wool alcohol ointment DAB	262037	1
Bombastus	Wool alcohol ointment DAB	263456	1
Bombastus	Wool alcohol ointment DAB	265139	1
Bombastus	Wool alcohol ointment DAB	265502	1
Bombastus	Wool alcohol ointment DAB	266204	1
Bombastus	Wool alcohol ointment DAB	267675	3
Noweda	Wool alcohol ointment DAB	267675	1
Gehe	Wool alcohol ointment DAB	268470	1
DAC 2011	Wool alcohol ointment DAB	268471	1
Bombastus	Wool alcohol ointment DAB	270640	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Wool alcohol ointment DAB	270642	1
Bombastus	Wool alcohol ointment DAB	270642	4
phönix 22.01.2013	Wool alcohol ointment DAB	270642	1
Noweda	Wool alcohol ointment DAB	270642	2
Bombastus	Wool alcohol ointment DAB	270780	2
Phönix	Wool alcohol ointment SR wit...	270781	1
Bombastus	Wool alcohol ointment SR wit...	270781	1
Bombastus	Wool alcohol ointment SR wit...	270784	2
Bombastus	Wool alcohol ointment DAB	272171	1
Sanacorp	Wool alcohol ointment DAB	274283	1
Defektur	Wool alcohol ointment DAB	2/13	1
Defektur	Wool alcohol ointment DAB	3/14	1
Apomix/Noweda	Wool alcohol ointment SR wit...	2A080713	1
Phönix Gotha	Wool alcohol ointment DAB	34554710	1
Caelo	Wool alcohol ointment DAB	4061204	1
Gehe 02.07.2013	Wool alcohol ointment DAB	4352610	1
Caelo	Wool alcohol ointment DAB	7121202	1
Apomix	Wool alcohol ointment SR wit...	6A110912	1
Defektur	Wool alcohol ointment DAB	6/14	1
Beiersdorf/ Phoenix	Wool alcohol ointment DAB	814Q1270	1
Beiersdorf/ Phoenix	Wool alcohol ointment DAB	814Q1370	1
Pharmachem/Phönix	Wool alcohol ointment DAB	90812	1
Caelo	Wool alcohol ointment DAB	11068001	1
Euro OTC	Wool alcohol ointment DAB	1108057-01	4
Euro OTC	Wool alcohol ointment DAB	110805701	1
Kehr	Wool alcohol ointment DAB	11184901	1
Caelo	Wool alcohol ointment DAB	11184901	1
Caelo	Wool alcohol ointment DAB	11184902	1
Caelo	Wool alcohol ointment DAB	11185001	1
Caesar & Loretz GmbH 40...	Wool alcohol ointment DAB	11185101	1
Caelo	Wool alcohol ointment DAB	11185201	1
Caelo	Wool alcohol ointment DAB	11232902	1
Kehr	Wool alcohol ointment DAB	11335001	1
Caelo	Wool alcohol ointment DAB	11335001	2
Caelo	Wool alcohol ointment DAB	11335003	1
Gehe	Wool alcohol ointment DAB	11335003	1
Phönix Gotha	Wool alcohol ointment DAB	11335003	2
Klenk	Wool alcohol ointment DAB	11335101	1
Anzag	Wool alcohol ointment DAB	11335101	1
Caelo	Wool alcohol ointment DAB	11335101	6
Kehr	Wool alcohol ointment DAB	1135001	1
Caelo	Wool alcohol ointment DAB	11335201	1
Caelo	Wool alcohol ointment DAB	12016202	3
Caelo	Wool alcohol ointment DAB	12016203	8
Caelo	Wool alcohol ointment DAB	12076601	3
Anzag Ek: 3,24EUR 11.09.12	Wool alcohol ointment DAB	12076601	1
Caelo	Wool alcohol ointment DAB	12076602	1
Caelo	Wool alcohol ointment DAB	12076701	2
Caelo	Wool alcohol ointment SR wit...	12078201	1
Noweda	Wool alcohol ointment SR wit...	12078201	1
Noweda	Wool alcohol ointment SR wit...	12078301	1
Euro OTC	Wool alcohol ointment DAB	1211018-01	1
Erika Apotheke	Wool alcohol ointment DAB	1215001	1
Caelo	Wool alcohol ointment DAB	12151801	2
Caelo	Wool alcohol ointment DAB	12151802	2
Caelo	Wool alcohol ointment DAB	12151901	8
Caelo	Wool alcohol ointment DAB	12151903	5
Bombastus	Wool alcohol ointment DAB	12151903	1
Phönix Gotha	Wool alcohol ointment DAB	12152001	2
Noweda	Wool alcohol ointment SR wit...	12244302	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol ointment DAB	12339901	4
Caelo	Wool alcohol ointment DAB	12340002	4
Caelo	Wool alcohol ointment DAB	1234002	1
Caelo	Wool alcohol ointment DAB	12340003	1
Caesar &Loretz GmbH	Wool alcohol ointment DAB	12340003	1
fragon	Wool alcohol ointment DAB	12345	2
Caelo	Wool alcohol ointment DAB	13066903	1
A.H. 05.06.13 ek: 3,33EUR	Wool alcohol ointment DAB	13067101	1
Caelo	Wool alcohol ointment DAB	13067101	2
Hustadt Apotheke	Wool alcohol ointment DAB	13067102	1
Caelo	Wool alcohol ointment DAB	13067201	11
Caelo	Wool alcohol ointment SR wit...	13077301	2
Caelo	Wool alcohol ointment SR wit...	13077302	1
Caelo	Wool alcohol ointment DAB	13093202	1
Alliance Healthcare	Wool alcohol ointment DAB	13114001	1
Caelo	Wool alcohol ointment DAB	13114103	4
Caelo	Wool alcohol ointment DAB	13126001	10
Gehe	Wool alcohol ointment DAB	13126001	1
Caelo	Wool alcohol ointment DAB	13126101	1
Caelo	Wool alcohol ointment DAB	13242401	1
Caelo	Wool alcohol ointment DAB	13355202	1
Caelo	Wool alcohol ointment DAB	1367101	1
Fagron	Wool alcohol ointment DAB	13D09-T04-003308	5
Fagron	Wool alcohol ointment DAB	13D09T04003308	2
Beiersdorf	Wool alcohol ointment DAB	14349010	1
Beiersdorf AG	Wool alcohol ointment DAB	14549310	1
Fagron	Wool alcohol ointment DAB	14819510	1
Caelo	Wool alcohol ointment DAB	20061204	1
Caelo	Wool alcohol ointment DAB	20271211-1	1
Beiersdorf	Wool alcohol ointment DAB	20271211-3	1
Beiersdorf AG	Wool alcohol ointment DAB	21111310	1
Beiersdorf / K	Wool alcohol ointment DAB	22142410	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	22412510	1
Caelo	Wool alcohol ointment DAB	22532810	1
Beiersdorf AG	Wool alcohol ointment DAB	22532810	1
Gehe	Wool alcohol ointment DAB	22532710	1
Beiersdorf/Otto	Wool alcohol ointment DAB	23033310	1
Gehe	Wool alcohol ointment DAB	23033310	1
Noweda	Wool alcohol ointment DAB	23223410	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	23223410	1
Phönix	Wool alcohol ointment DAB	23733710	1
Caelo	Wool alcohol ointment DAB	24655410	1
Caelo	Wool alcohol ointment DAB	25051306	1
Bombastus	Wool alcohol ointment SR wit...	270776	1
Bombastus	Wool alcohol ointment SR wit...	274517	2
Bombastus	Wool alcohol ointment SR wit...	274518	3
Bombastus	Wool alcohol ointment SR wit...	274521	1
Bombastus	Wool alcohol ointment DAB	275999	5
Phönix	Wool alcohol ointment DAB	275999	1
Bombastus	Wool alcohol ointment SR wit...	276135	1
Bombastus	Wool alcohol ointment SR wit...	277023	2
Gehe	Wool alcohol ointment DAB	277023	1
Bombastus	Wool alcohol ointment DAB	277190	1
Bombastus	Wool alcohol ointment DAB	277191	1
Gehe	Wool alcohol ointment DAB	277191	2
Bombastus	Wool alcohol ointment DAB	277192	3
Bombastus	Wool alcohol ointment DAB	278684	3
Bombastus	Wool alcohol ointment SR wit...	279104	1
Bombastus	Wool alcohol ointment SR wit...	279105	1
Bombastus	Wool alcohol ointment DAB	279934	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol ointment DAB	280057	1
Bombastus	Wool alcohol ointment DAB	280058	1
Bombastus	Wool alcohol ointment SR wit...	280060	1
Gehe	Wool alcohol ointment DAB	280061	2
Bombastus	Wool alcohol ointment SR wit...	280061	2
Fagron	Wool alcohol ointment DAB	281690	1
Bombastus	Wool alcohol ointment DAB	283175	1
Phoenix	Wool alcohol ointment DAB	283619	1
Caelo	Wool alcohol ointment DAB	30111302	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	30711310	1
Beiersdorf AG	Wool alcohol ointment DAB	30721310	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	30731210	1
Phönix	Wool alcohol ointment DAB	30821610	1
Caelo	Wool alcohol ointment DAB	30831410	1
AHD, 05.04.2013	Wool alcohol ointment DAB	30831410	1
Phoenix	Wool alcohol ointment DAB	30831410	1
Beiersdorf/Phönix	Wool alcohol ointment DAB	30841410	1
Beiersdorf	Wool alcohol ointment DAB	30851510	1
Caelo	Wool alcohol ointment DAB	31251910	1
Phönix Gotha	Wool alcohol ointment DAB	31231810	1
Beiersdorf AG	Wool alcohol ointment DAB	31251910	1
Beiersdorf	Wool alcohol ointment DAB	31322010	1
Großhandel	Wool alcohol ointment DAB	31322010	1
Beiersdorf/ Noweda	Wool alcohol ointment DAB	31322010	1
Gehe	Wool alcohol ointment DAB	31322010	1
BDF Beiersdorf AG	Wool alcohol ointment DAB	31322010	1
Phoenix	Wool alcohol ointment DAB	32344210	1
Beiersdorf/Gehe	Wool alcohol ointment DAB	32352710	1
Anzag	Wool alcohol ointment DAB	3312E-01370	1
Beiersdorf AG	Wool alcohol ointment DAB	33323210	1
Beiersdorf / Noweda	Wool alcohol ointment DAB	33343310	1
BDF/ Noweda	Wool alcohol ointment DAB	33343310	1
Beiersdorf/Phoenix	Wool alcohol ointment DAB	33723610	1
eucerinum	Wool alcohol ointment DAB	33753710	1
Beiersdorf AG	Wool alcohol ointment DAB	34344210	1
Noweda	Wool alcohol ointment DAB	34414610	1
Beiersdorf AG/Noweda	Wool alcohol ointment DAB	34414610	1
eucerinum	Wool alcohol ointment DAB	34554710	1
Defektur	Wool alcohol ointment DAB	3812	1
Defektur	Wool alcohol ointment DAB	40/13	1
Caelo	Wool alcohol ointment DAB	44/AKD9	1
Defektur	Wool alcohol ointment DAB	47/13	1
Phönix/beiersdorf	Wool alcohol ointment DAB	50955510	1
Defektur	Wool alcohol ointment DAB	53/13	1
Caelo	Wool alcohol ointment DAB	U-240513-3	1

- 81 792 spectra from 499 *Apo-Ident* customers from a total of 24 501 batches from a further 322 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Paraffin classification* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose,

all relevant spectra of the various substances were compared with *Paraffin classification* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	2489	0	33 055
Type B	0	5599	10	88 067
Type C	0	9427	133	81 750

The substance/substance group *Paraffin classification* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9702 %)	100.0000 % (> 99.7589 %)
Type B	100.0000 % (> 99.9764 %)	99.8217 % (> 99.7682 %)
Type C	100.0000 % (> 99.6595 %)	98.6088 % (> 98.5774 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30874	31702	3.32	21.45
30877	31358	5.25	24.85
30932	30932	0.00	17.52
30933	30933	0.00	15.96
31009	31941	3.11	27.54
31010	31941	4.69	28.44
31168	31446	3.04	23.10
31174	31174	0.00	26.73
31228	31228	0.00	30.92
31272	31272	0.00	30.69
31292	31446	3.53	23.62
31351	31351	0.00	19.26
31358	31358	0.00	25.99
31437	31916	2.16	26.63
31441	32017	8.15	31.55
31445	31610	3.92	27.74
31446	31446	0.00	22.28
31456	31456	0.00	34.59
31468	31916	3.97	26.46
31473	31473	0.00	25.43

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31522	31522	0.00	30.23
31541	31736	9.22	28.55
31567	31456	8.42	35.90
31570	31570	0.00	19.78
31587	31702	3.26	23.15
31590	31736	7.98	30.04
31601	32495	3.34	26.84
31602	31602	0.00	35.18
31610	31610	0.00	27.34
31611	31909	6.20	37.60
31612	31456	6.23	36.03
31702	31702	0.00	21.02
31736	31736	0.00	28.55
31740	31740	0.00	30.45
31744	31744	0.00	28.95
31807	31602	4.79	35.87
31812	31812	0.00	29.32
31816	31816	0.00	15.85
31909	31909	0.00	36.41
31910	31910	0.00	38.61
31916	31916	0.00	24.99
31917	31917	0.00	27.33
31920	31920	0.00	15.99
31921	31921	0.00	30.61
31940	31940	0.00	22.25
31941	31941	0.00	23.75
31949	31744	1.98	27.39
31954	31602	2.67	34.80
31999	31999	0.00	26.79
32017	32017	0.00	36.61
32279	32279	0.00	17.04
32495	32495	0.00	30.54
32517	32517	0.00	31.55
32704	32704	0.00	30.27
32900	32900	0.00	29.14
32946	31174	6.18	24.80
32947	31999	8.24	25.74
33133	33133	0.00	25.31
33274	32704	3.67	27.75
33396	33396	0.00	16.41
33729	31702	3.76	22.14

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Pasta Cordes®**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30633-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Pasta Cordes®

Special notes

When selecting the *Pasta Cordes®* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Pasta Cordes®*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Ichthyol	Pasta Cordes®	13A056	31084	40	1406053
Ichthyol	Pasta Cordes®	16A063	33285	40	20160427

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Pasta Cordes*[®]. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 3 reference samples of the substance/substance group *Pasta Cordes*[®].
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Ichthyol	Pasta Cordes [®]	12A008	30633	40
Ichthyol	Pasta Cordes [®]	13A056	31673	60
Ichthyol	Pasta Cordes [®]	17A173	34046	40

- 93 734 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 175 spectra from 54 *Apo-Ident* customers from 29 batches from the substance/substance group *Pasta Cordes*[®].
- These include spectra of independent samples from 27 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
ichthyol/Gehe	Pasta Cordes [®]	12A182	2
Ichthyol /Gehe	Pasta Cordes [®]	12A182	1
Ichthyol/Sanacorp	Pasta Cordes [®]	12A182	1
Ichthyol /Noweda	Pasta Cordes [®]	13A056	3
Ichthyol	Pasta Cordes [®]	13A056	1
Ichthyol Gesellschaft	Pasta Cordes [®]	13A056	2
Noweda	Pasta Cordes [®]	13A056	1
Ichthyol-Gesellschaft	Pasta Cordes [®]	13A056	1
Caelo	Pasta Cordes [®]	13A098	3
Caelo	Pasta Cordes [®]	13A056	4
Ichthyol/Fiebig	Pasta Cordes [®]	13A056	1
Ichthyol/Gehe	Pasta Cordes [®]	13A098	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	13A098	1
Ichthyol	Pasta Cordes [®]	13A098	3
Ichthyol/Ebert	Pasta Cordes [®]	13A098	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	13A146	3
ichthyol/Gehe	Pasta Cordes [®]	13A146	1
Ichthyol/Phönix	Pasta Cordes [®]	14A025	2
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	14A025	2
Cordes/Gehe	Pasta Cordes [®]	14A025	1
Ichthyol	Pasta Cordes [®]	14A025	1
Ichthyol/Phoenix	Pasta Cordes [®]	14A025	1
Ichthyol/Gehe	Pasta Cordes [®]	14A124	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	14A124	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	14A124	1
Ichthyol/Phönix	Pasta Cordes [®]	14A124	1
Ichthyol	Pasta Cordes [®]	14A124	2
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	14a124	1
Ichthyol Gesellschaft C...	Pasta Cordes [®]	14A124	1
Sanacorp / Ichthyolgese...	Pasta Cordes [®]	14A124	1
Ichthyol-Gesellschaft	Pasta Cordes [®]	14A124	1
Ichthyol/Noweda	Pasta Cordes [®]	14A142	2
Ichthyol/ Phönix	Pasta Cordes [®]	14A142	1
Ichthyol/Phönix	Pasta Cordes [®]	14A142	7
Ichthyol/Fiebig	Pasta Cordes [®]	14A142	2
Ichthyol	Pasta Cordes [®]	14A142	1
Fischer/Gehe	Pasta Cordes [®]	14A142	2
Ichthyol Gesellschaft	Pasta Cordes [®]	14A142	1
Ichthyol / AHD	Pasta Cordes [®]	14A142	1
Ichthyol / Noweda	Pasta Cordes [®]	14A142	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	14A199	4
Ichthyol/Gehe	Pasta Cordes [®]	14A199	14
Ichthyol	Pasta Cordes [®]	14A199	1
Phönix/Ichthyol	Pasta Cordes [®]	15A121	2
Ichthyol/ Noweda	Pasta Cordes [®]	14A025	1
ICHTHYOL-GESELLSCHAFT	Pasta Cordes [®]	15A121	1
Ichthyol Gesellschaft/K...	Pasta Cordes [®]	15A121	1
Noweda-Ichthyol	Pasta Cordes [®]	15A121	1
Ichthyol/Phoenix	Pasta Cordes [®]	15A121	1
Ichthyol/Sanacorp	Pasta Cordes [®]	15A121	1
Ichthyol Gesellschaft/N...	Pasta Cordes [®]	15A121	1
ICHTYOL GESELLSCHAFT Co...	Pasta Cordes [®]	15A121	1
Ichthyol	Pasta Cordes [®]	15A121	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe/Ichthyol	Pasta Cordes [®]	15A173	1
Ichthyol/Sanacorp	Pasta Cordes [®]	15A173	2
Ichthyol/Phönix WE:2804...	Pasta Cordes [®]	15A121	1
Ichthyol Cordes/Phoenix	Pasta Cordes [®]	15A173	1
Ichthyol/Noweda	Pasta Cordes [®]	15A173	1
Ichthyol-Gesellschaft	Pasta Cordes [®]	15A121	1
Ichthyol/ Krieger	Pasta Cordes [®]	15A173	1
Cordes / Phönix	Pasta Cordes [®]	15A173	1
Ichthyol/Phoenix	Pasta Cordes [®]	15A173	1
WE 28.11.16 Cordes/Phoe...	Pasta Cordes [®]	15A173	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	15A173	2
Ichthyol/OG	Pasta Cordes [®]	15A173	1
Ichthyol	Pasta Cordes [®]	15A173	1
Caelo	Pasta Cordes [®]	160822PC	2
Ichthyol/Phoenix	Pasta Cordes [®]	16A063	1
Caelo	Pasta Cordes [®]	16A063	2
Ichthyol Gesellschaft	Pasta Cordes [®]	16A063	1
Ichthyol/Noweda	Pasta Cordes [®]	16A063	1
Ichthyol/Gehe	Pasta Cordes [®]	16A064	21
Ichthyol/AHD	Pasta Cordes [®]	16A063	1
Ichthyol	Pasta Cordes [®]	16A064	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	16A064	2
Ichthyol/Noweda	Pasta Cordes [®]	16A098	1
Ichthyol Gesellschaft C...	Pasta Cordes [®]	16A098	1
Ichthyol	Pasta Cordes [®]	17A173	1
Ichthyol-Gesellschaft	Pasta Cordes [®]	170102PC/15A173	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	17A173	1
Ichthyol-Gesellschaft C...	Pasta Cordes [®]	17A021	1
Ichthyol/Phoenix	Pasta Cordes [®]	17A021	1
Ichthyol Gesellschaft /...	Pasta Cordes [®]	180308PC/16A064	1
ICHTHYOL	Pasta Cordes [®]	10A187	1
Ichthyol Ges	Pasta Cordes [®]	11A069	1
Caelo	Pasta Cordes [®]	11A206	2
Ichthyol/Gehe	Pasta Cordes [®]	11A206	2
Ichthyol	Pasta Cordes [®]	12A008	1
Ichthyol/Gehe	Pasta Cordes [®]	12A043	3
Ichthyol Ges	Pasta Cordes [®]	12A138	1
Phoenix	Pasta Cordes [®]	12A043	1
Ichthyol/Sanacorp	Pasta Cordes [®]	12A138	1
Phönix	Pasta Cordes [®]	12A182	1
Ichthyol	Pasta Cordes [®]	12a043	1
Ichthyol	Pasta Cordes [®]	12a182	2
Ichthyol/Phoenix	Pasta Cordes [®]	12A182	1
Ichthyol /Noweda	Pasta Cordes [®]	12A182	1
Caelo	Pasta Cordes [®]	12A182	2
Ichthyol / Gehe	Pasta Cordes [®]	12A182	1
Caelo	Pasta Cordes [®]	12AD43	1
ichthyol gesellschaft	Pasta Cordes [®]	13A098	1
Cordes/AHD	Pasta Cordes [®]	13A056	1

- 91 177 spectra from 501 *Apo-Ident* customers from a total of 26 980 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Pasta Cordes*[®] can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Pasta Cordes*[®] and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	126	14	93 734
Type C	0	142	33	91 152

The substance/substance group *Pasta Cordes*[®] can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	90.0000 % (> 87.8571 %)
Type C	100.0000 % (> 99.6586 %)	81.1429 % (> 79.4286 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31084	31084	0.00	35.31
33285	31084	4.85	35.88

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Permethrin 25% solution
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31177-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Permethrin 25% solution

Special notes

When selecting the *Permethrin 25% solution* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Permethrin 25% solution*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Infectopharm	Permethrin 25% s...	S031301.1	31177	40	1405128
Infectopharm	Permethrin 25% s...	S061401.1	31563	40	1508108
Infectopharm	Permethrin 25% s...	S021601.1	32262	40	20160310
Infectopharm	Permethrin 25% s...	S021601.1	32263	40	20160310
Infectopharm	Permethrin 25% s...	S071602.1	33247	40	20160728

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Infectopharm	Permethrin 25% s...	S071602.1	33301	40	20160728
Infectopharm	Permethrin 25% s...	S021801.1	34120	40	20180301
Infectopharm	Permethrin 25% s...	S021801.1	34121	40	20180301

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 320 spectra of 8 reference samples from the substance/substance group *Permethrin 25% solution*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 224 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 300 spectra from 8 reference samples of the substance/substance group *Permethrin 25% solution*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Infectopharm	Permethrin 25% solution	S061101.1	30583	40
Infectopharm	Permethrin 25% solution	S021201.1	30584	40
Infectopharm	Permethrin 25% solution	S061401.1	31563 [†]	20
Infectopharm	Permethrin 25% solution	S071602.1	33246	40
Infectopharm	Permethrin 25% solution	S071602.1	33302	40
Infectopharm	Permethrin 25% solution	S071602.1	33328	40
Infectopharm	Permethrin 25% solution	S041701.2	33940	40
Infectopharm	Permethrin 25% solution	S041701.2	33941	40

- 93 574 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 173 spectra from 75 *Apo-Ident* customers from 52 batches from the substance/substance group *Permethrin 25% solution*.
- These include spectra of independent samples from 30 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
InfectoPharm	Permethrin 25% solution	1072015A	1
InfectoPharm	Permethrin 25% solution	2091402	1
Infectopharm/noweda	Permethrin 25% solution		1
InfectoPharm	Permethrin 25% solution	4071506	1
InfectoPharm	Permethrin 25% solution	5011501	1
InfectoPharm	Permethrin 25% solution	9101406	1
InfectoPharm	Permethrin 25% solution	11091403	1
Caelo	Permethrin 25% solution		6
InfectoPharm	Permethrin 25% solution	14041401	1
Biomo/Noweda	Permethrin 25% solution	150901	1
Biomo/Noweda	Permethrin 25% solution	160202	2
InfectoPharm	Permethrin 25% solution	25111401	1
InfectoPharm	Permethrin 25% solution	28011504	1
InfectoPharm	Permethrin 25% solution	20081401	1
InfectoPharm	Permethrin 25% solution	31071402	1
infectopharm	Permethrin 25% solution	23071414	1
InfectoPharm	Permethrin 25% solution	17121403	1
Infectopharm/Phoenix	Permethrin 25% solution		4
Phoenix	Permethrin 25% solution		3
InfectoPharm	Permethrin 25% solution		9
Infectopharm	Permethrin 25% solution		27
Fagron	Permethrin 25% solution		1
InfectoPharm / Phönix	Permethrin 25% solution		2
AHD / Infectopharm	Permethrin 25% solution		1
Infectopharm/ Noweda	Permethrin 25% solution		7
infectopharm	Permethrin 25% solution		2
infectopharm / Sanacorp	Permethrin 25% solution		1
InfectoPharm	Permethrin 25% solution	S0216011	1
AHD / Infectopharm	Permethrin 25% solution	S0216011	1
Infectopharm Arzneimittel	Permethrin 25% solution	S0216011	1
Infectopharm/Phönix	Permethrin 25% solution		3
	Permethrin 25% solution		1
Infectopharm	Permethrin 25% solution	S0614011	1
Euro OTC	Permethrin 25% solution		1
INFECTOPHARM	Permethrin 25% solution		2

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Infectopharm/Gehe	Permethrin 25% solution	S0313011	1
Noweda/InfectoPharm	Permethrin 25% solution		1
Gehe Infectopharma	Permethrin 25% solution		1
Noweda/ Infectopharm	Permethrin 25% solution		1
InfectoPharm/AHD	Permethrin 25% solution		1
InfectoPharm/Kehr	Permethrin 25% solution		1
Bombastus	Permethrin 25% solution		1
Infectopharm Arzneimittel	Permethrin 25% solution		2
Infectopharm/Anzag	Permethrin 25% solution		1
Infectopharm/Gehe	Permethrin 25% solution		1
INFECTOPHARM/Sanacorp	Permethrin 25% solution		1
Infecto Pharm / Phönix	Permethrin 25% solution	S072018	1
Phoenix/InfectoPharm	Permethrin 25% solution		5
Infectopharm/AHD	Permethrin 25% solution		6
Infectopharm / 20	Permethrin 25% solution		1
InfectoPharm/Phönix	Permethrin 25% solution		1
InfectoPharm/Gehe	Permethrin 25% solution		2
Infectopharm / Alliance...	Permethrin 25% solution		1
Infectopharm/Ebert	Permethrin 25% solution		1
Infectopharm/Sanacorp	Permethrin 25% solution		4
Infectopharm/ Sanacorp	Permethrin 25% solution		3
Infectopharm/Noweda	Permethrin 25% solution		2
InfectoPharma / Noweda	Permethrin 25% solution		1
Infectopharm / Kehr	Permethrin 25% solution		2
Dermapharm	Permethrin 25% solution		2
Ichthyol/Phönix	Permethrin 25% solution		2
Infecto Pharm	Permethrin 25% solution		2
Infect Pharm	Permethrin 25% solution		2
Infectopharm/Sana	Permethrin 25% solution		1
Phönix/Infectopharm	Permethrin 25% solution		1
Infectopharm	Permethrin 25% solution	S1214021	1
InfectoPharm	Permethrin 25% solution	6011407	1
InfectoPharm	Permethrin 25% solution	8011424	1
InfectoPharm	Permethrin 25% solution	11061302	1
InfectoPharm	Permethrin 25% solution	21121202	1
InfectoPharm	Permethrin 25% solution	23111205	1
InfectoPharm	Permethrin 25% solution	25021307	1
Anzag (Hersteller: Infe...	Permethrin 25% solution		1
InfectoPharm	Permethrin 25% solution	2510210212011	1
Infectopharm/ Phoenix	Permethrin 25% solution		1
Infectopharm, von Humbo...	Permethrin 25% solution		3
Infectopharm/ Now	Permethrin 25% solution		2
InfectoPharm/ G	Permethrin 25% solution		1
Infecto/Phönix	Permethrin 25% solution		1
Infectopharm	Permethrin 25% solution	S0313011	1
infecto Pharma	Permethrin 25% solution		4
Löwen City-Apotheke Bar...	Permethrin 25% solution		1
Noweda	Permethrin 25% solution		2
Gehe	Permethrin 25% solution		1
Infectopharm	Permethrin 25% solution	S081202	1
infectopharm/Jenne	Permethrin 25% solution		1
Infectopharm/Phönix	Permethrin 25% solution	S0313011	1
infecto pharm	Permethrin 25% solution	s121402	1

- 91 179 spectra from 501 *Apo-Ident* customers from a total of 26 958 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Permethrin 25% solution* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Permethrin 25% solution* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	320	0	35 224
Type B	0	290	10	93 574
Type C	0	165	8	91 179

The substance/substance group *Permethrin 25% solution* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.1250 %)
Type B	100.0000 % (> 99.9761 %)	96.6667 % (> 95.6667 %)
Type C	100.0000 % (> 99.6586 %)	95.3757 % (> 93.6416 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31177	31177	0.00	123.96
31563	31563	0.00	135.97
32262	31563	7.28	164.56
32263	31563	4.79	159.34
33247	31177	8.45	136.36
33301	31177	18.70	140.13
34120	31177	19.53	147.64
34121	31177	14.18	150.07

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Petroleum ether**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30899-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Petroleum ether; Aether petrolei; Benzinum (petrolei)

Special notes

When selecting the *Petroleum ether* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Petroleum ether*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Petroleum ether	12279604	30899	60	1503677
Otto Fischar	Petroleum ether	6323042	31253	40	1505436
Caelo	Petroleum ether	13141304	31283	40	1412255

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 140 spectra of 3 reference samples from the substance/substance group *Petroleum ether*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 404 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Petroleum ether*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Petroleum ether	91883269	30102	40
Hedinger	Petroleum ether	028175B	32955	40
Caelo	Petroleum ether	171715	33329	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1159 spectra from 257 *Apo-Ident* customers from 254 batches from the substance/substance group *Petroleum ether*.
- These include spectra of independent samples from 247 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Hedinger/Gehe	Petroleum ether	28170B	2
Hedinger	Petroleum ether	28171B	25
Hedinger vdL	Petroleum ether	28171B	2
Hedinger/Anzag	Petroleum ether	28171B	2
Hedinger	Petroleum ether	28172b	1
Hedinger/AHD	Petroleum ether	28174B	1
Hedinger/Phönix	Petroleum ether	28174B	2
Hedinger	Petroleum ether	28174B	20
Hedinger/Phoenix	Petroleum ether	28175B	1
Hedinger/Noweda	Petroleum ether	28175B	4
HEDINGER / NOWEDA	Petroleum ether	28175B	1
Otto Fischar GmbH	Petroleum ether	6324063	1
Fischar/Spangro	Petroleum ether	6326026	1
Caelo	Petroleum ether	6326031	1
Hedinger	Petroleum ether	28168B	19
Fischar	Petroleum ether	6323097	1
Hedinger	Petroleum ether	28169B	13
Fischar/Phoenix	Petroleum ether	6326047	1
Hedinger/AHD	Petroleum ether	28167B	3
Fischar/Gehe	Petroleum ether	6323093	1
Fischar GmbH	Petroleum ether	6323066	1
	Petroleum ether		1
Hedinger/ Phönix	Petroleum ether	(028)028	1
Hedinger	Petroleum ether	281658	1
Hedinger	Petroleum ether	2071202	1
Noweda	Petroleum ether	281658	1
Hedinger	Petroleum ether	28165B	6
Caelo	Petroleum ether	281668	1
Hedinger, Heiligenwiese...	Petroleum ether	28166B	2
Hedinger	Petroleum ether	28166B	9
Hedinger	Petroleum ether	281678	1
Hedinger/Noweda	Petroleum ether	281678	1
Sanacorp	Petroleum ether	28166B	1
Hedinger/Otto	Petroleum ether	281678	1
Klenk	Petroleum ether	28167b	1
Hedinger	Petroleum ether	28167B	18
Noweda/Hedinger	Petroleum ether	28166B	1
Hedinger / Gehe	Petroleum ether	28167B	1
Caelo	Petroleum ether	28166B	1
Hedinger Phönix	Petroleum ether	28167B	1
Hedinger/Stumpf	Petroleum ether	28167B	1
Hedinger/Anzag	Petroleum ether	28167B	2
Hedinger/Anzag 20.11.2012	Petroleum ether	28167B	2
Hedinger/Gehe	Petroleum ether	28167B	2
Gehe/ Hedinger	Petroleum ether	28167B	1
heidinger/gehe	Petroleum ether	28167B	1
Anzag/Hedinger	Petroleum ether	28167B	1
Hediniger	Petroleum ether	28167B	1
HEDINGER	Petroleum ether	281688	2
hedinger/phönix	Petroleum ether	28167b	1
Hedinger	Petroleum ether	281688	3
hedinger	Petroleum ether	28168B	1
Caelo	Petroleum ether	28168B	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Gehe	Petroleum ether	28168B	2
Hedinger/ jenne	Petroleum ether	28168B	3
Hedinger/Sanacorp	Petroleum ether	28168B	3
Hedinger/Phoenix	Petroleum ether	28168B	1
Caelo	Petroleum ether	281698	1
Noweda	Petroleum ether	281698	1
hedinger/noweda	Petroleum ether	281698	1
Hedinger	Petroleum ether	281698	2
Hedinger/AHD	Petroleum ether	28169B	1
hedinger/anzag	Petroleum ether	28169B	1
Hedinger / Gehe	Petroleum ether	28169B	3
Spangr/Hedinger	Petroleum ether	28169B	1
Caelo	Petroleum ether	28169B	9
hedinger/noweda	Petroleum ether	28169B	1
Hedinger / NOW	Petroleum ether	28170B	1
Hedinger / Sanacorp	Petroleum ether	28169B	2
Sanacorp	Petroleum ether	28169B	1
Hedinger	Petroleum ether	28170B	14
Hedinger , Stuttgart	Petroleum ether	28170B	1
Hedinger / Gehe	Petroleum ether	28170B	5
Hedinger/Sanacorp	Petroleum ether	28170B	1
Hedinger/Jenne	Petroleum ether	28170B	4
Hedinger/V.d.L.	Petroleum ether	28170B	1
HEDINGER / 20	Petroleum ether	28170B	1
Hedinger/Noweda	Petroleum ether	28170B	2
Anzag 16.10.14	Petroleum ether	28171B	1
Sanacorp	Petroleum ether	28171B	5
Barthel	Petroleum ether	28171B	3
Hedinger/GEHE	Petroleum ether	28171B	1
Hedinger/Phönix	Petroleum ether	28171B	5
Hedinger/ Noweda	Petroleum ether	28171B	1
Hedinger / Gehe	Petroleum ether	28171B	1
Hedinger/Gehe	Petroleum ether	28171B	6
Caelo	Petroleum ether	28171B	4
Gehe - Hedinger	Petroleum ether	28171B	1
hedinger/gehe	Petroleum ether	28171B	1
Spangro/Hedinger	Petroleum ether	28171B	1
Hedinger/Noweda	Petroleum ether	28171B	2
A. Hedinger / Noweda	Petroleum ether	281728	1
Hedinger	Petroleum ether	28172B	37
Caelo	Petroleum ether	28172B	2
Hedinger / Phönix	Petroleum ether	28172B	1
Hedinger/Sanacorp	Petroleum ether	28172B	3
Hedinger / Noweda	Petroleum ether	28172B	2
Hedinger/Sanacorp	Petroleum ether	28172b	1
Hedinger/ Noweda	Petroleum ether	28172B	1
Hedinger/Geilenkirchen	Petroleum ether	28172B	2
phoenix/hedinger	Petroleum ether	28172B	1
heidinger/gehe	Petroleum ether	28172B	2
Hedinger/Gehe	Petroleum ether	28172B	6
Hedinger/AllianceHC	Petroleum ether	28172B	1
Hedinger	Petroleum ether	281738	4
Caelo	Petroleum ether	28173B	4
Hedinger	Petroleum ether	28173B	15
Hedinger/Gehe	Petroleum ether	28173B	5
Gehe	Petroleum ether	28173b	1
Hedinger/AHD	Petroleum ether	28173B	1
Hedinger/Sanacorp	Petroleum ether	28173b	1
Hedinger/Noweda	Petroleum ether	28173B	3
Hedinger	Petroleum ether	281748	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/AllianceHC	Petroleum ether	28173B	1
Hedinger/Alliance Healt...	Petroleum ether	28174B	1
Caelo	Petroleum ether	28174B	2
Hedinger Phönix	Petroleum ether	28174B	1
Hedinger/Noweda	Petroleum ether	28174B	6
Hedinger/Sanacorp	Petroleum ether	28174B	1
Hedinger/Gehe	Petroleum ether	28174B	1
Hedinger	Petroleum ether	28174b	2
heidinger/anzag	Petroleum ether	28174B	1
Hedinger/ Gehe	Petroleum ether	28174B	1
Hedinger	Petroleum ether	28175B	14
Hedinger/Sanacorp	Petroleum ether	28175B	2
Caelo	Petroleum ether	28175B	6
Hedinger/Gehe	Petroleum ether	28175B	4
heidinger/gehe	Petroleum ether	28175B	1
Calo/AHD	Petroleum ether	28175B	1
L 20/Hedinger	Petroleum ether	28175B	1
Gehe/Hedinger	Petroleum ether	28175B	1
Hedinger GmbH	Petroleum ether	28175B	1
Sanacorp	Petroleum ether	28175B	1
Noweda/hedinger	Petroleum ether	28175B	1
Hedinger/Geilenkirchen	Petroleum ether	28176B	1
Hedinger	Petroleum ether	28176B	5
Hedinger	Petroleum ether	281768	3
heidinger	Petroleum ether	28176B	1
GEHE	Petroleum ether	28176B	1
Hedinger/Sanacorp	Petroleum ether	28176B	1
Caelo	Petroleum ether	28176B	5
Hedinger/Kehr	Petroleum ether	28176B	1
Hedinger/noweda	Petroleum ether	28176B	1
Hedinger/Phoenix	Petroleum ether	28177B	1
Hedinger	Petroleum ether	2817SB	1
Hedinger	Petroleum ether	3071304	1
Hedinger	Petroleum ether	3071401	1
Caelo	Petroleum ether	4031404	1
Otto Fischer GmbH	Petroleum ether	15031402	1
Caelo	Petroleum ether	11357104	1
otto fischer	Petroleum ether	15081602	1
Hedinger/GEHE	Petroleum ether	1914M1170	1
Jenne	Petroleum ether	15045203	1
Fischar	Petroleum ether	24061501	1
Hedinger	Petroleum ether	2411Q-01170	1
Hedinger/ Noweda	Petroleum ether	281698	1
Hedinger	Petroleum ether	28011506	1
Hedinger	Petroleum ether	25021302	1
Noweda/Hedinger	Petroleum ether	281748	1
Hedinger	Petroleum ether	30101310	1
Hedinger	Petroleum ether	31101203	1
Hedinger GmbH/Noweda	Petroleum ether	281738	1
Hedinger	Petroleum ether	29041501	1
Fischar / Phönix	Petroleum ether	3413J-01170	1
Fischar / Phönix	Petroleum ether	3412I-01170	1
Fischar/Gehe	Petroleum ether	3623111	1
Hedinger / Phoenix	Petroleum ether	3966/V02/17	1
Fischar	Petroleum ether	412062-2-P	1
Fischar	Petroleum ether	55780	1
fischar	Petroleum ether	616331022	1
fischar	Petroleum ether	6320031	1
Phönix	Petroleum ether	6321111	1
Otto Fischar GmbH	Petroleum ether	6321091	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischar	Petroleum ether	6321111	1
Fischar	Petroleum ether	6322041	1
Fischer Saarbrücken/GEHE	Petroleum ether	6322041	1
Phönix/Fischar	Petroleum ether	6322044	6
Otto Fischar GmbH	Petroleum ether	6322045	1
Fischar/Gehe	Petroleum ether	6322053	1
Fischar	Petroleum ether	6322051	1
Fischar	Petroleum ether	6322061	1
Gehe	Petroleum ether	6322061	1
Fischar/Gehe	Petroleum ether	6322061	1
Phoenix	Petroleum ether	6322074	1
Fischar	Petroleum ether	6322074	2
Fischar	Petroleum ether	6322079	1
Fischar	Petroleum ether	6322076	1
Fischar	Petroleum ether	632207A	3
Fischar / Anzag	Petroleum ether	6322095	1
Gehe	Petroleum ether	6322095	1
Sanacorp	Petroleum ether	6322096	1
Caelo	Petroleum ether	6322097	1
Anzag/Fischar	Petroleum ether	6322098	1
Phönix	Petroleum ether	6322101	9
Phönix/Fischar	Petroleum ether	6322102	16
Phönix/Fischar	Petroleum ether	6322103	1
Fischer/Phoenix	Petroleum ether	6322104	1
Phönix/Fischar	Petroleum ether	6322112	6
Fischar / Noweda	Petroleum ether	6322112	1
Großhandel	Petroleum ether	6322113	1
Fischar	Petroleum ether	6322113	1
Phönix/Fischar	Petroleum ether	6322113	5
Fischar	Petroleum ether	6322114	1
Fischar/Sanacorp	Petroleum ether	6322116	1
Otto Fischar	Petroleum ether	6322116	1
Fischar/Gehe	Petroleum ether	6322117	1
gehe, 14.01.2013	Petroleum ether	6322117	1
Phönix	Petroleum ether	6322118	6
GEHE	Petroleum ether	6322119	1
Fischar	Petroleum ether	6322119	1
Fagron	Petroleum ether	6323014	3
Phönix/ Fischar	Petroleum ether	6323015	12
Fischar/Phoenix	Petroleum ether	6323016	1
Fischar	Petroleum ether	6323016	1
Fischar	Petroleum ether	6323017	1
Fischar/Noweda	Petroleum ether	6323017	1
Anzag/Fischar	Petroleum ether	6323017	1
Fischar GmbH	Petroleum ether	6323022	1
Großhandel	Petroleum ether	6323032	1
Caelo	Petroleum ether	6323033	2
Fischar	Petroleum ether	6323042	11
Fischer Saarbrücken/GEHE	Petroleum ether	6323043	1
Fischar / Sanacorp	Petroleum ether	6323062	1
Fischar	Petroleum ether	6323063	3
Fischar/Anzag	Petroleum ether	6323062	1
Fischar/Phönix	Petroleum ether	6323063	1
Otto Fischar	Petroleum ether	6323063	1
fischar/phoenix	Petroleum ether	6323064	1
Fischar	Petroleum ether	6323067	6
Fischer	Petroleum ether	6323064	1
Fischar/ Now	Petroleum ether	6323064	1
Fischar/Noweda	Petroleum ether	6323065	1
Fischar	Petroleum ether	6323066	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
otto fischar/Sana	Petroleum ether	6323065	1
Caelo	Petroleum ether	6323068	17
Sanacorp	Petroleum ether	6323068	1
Fischar	Petroleum ether	6323068	1
Otto Fischer	Petroleum ether	6323069	1
Fischar/Phoenix	Petroleum ether	6323069	1
Otto Fischar Phönix	Petroleum ether	6323069	1
Fischar / Noweda	Petroleum ether	6323069	1
Fischar	Petroleum ether	6323091	6
Fischar/Noweda	Petroleum ether	6323069	1
Ebert/Fischar	Petroleum ether	6323093	1
Fischar	Petroleum ether	6323093	7
FISCHER GMBH / GEHE	Petroleum ether	6323094	1
Otto Fischar GmbH	Petroleum ether	6323094	1
Caelo	Petroleum ether	6323095	1
Fischar/Anzag	Petroleum ether	6323099	1
Fischar/Anzag	Petroleum ether	6323095	1
Fischar	Petroleum ether	6323099	1
Otto Fischar Phönix	Petroleum ether	6323111	1
Fischar	Petroleum ether	6323111	1
Fischer/Sanacorp	Petroleum ether	6323099	1
Großhandel	Petroleum ether	6323111	1
Fischar	Petroleum ether	6323112	5
Fischar/Nowenda	Petroleum ether	6323112	1
Fischar	Petroleum ether	6324013	1
Fischar	Petroleum ether	6324016	1
Caelo	Petroleum ether	6324017	29
Fischar/Gehe	Petroleum ether	6323113	1
Fischar Phönix	Petroleum ether	6324018	1
Caelo	Petroleum ether	6324018	5
Fischar/Phoenix	Petroleum ether	6324018	1
Fischar	Petroleum ether	632401A	1
Alliance Healthcare/Fis...	Petroleum ether	6324025	1
Fischar/Sanacorp	Petroleum ether	6324033	1
Fischar/ Phönix	Petroleum ether	6324034	13
Fischar	Petroleum ether	6324034	2
Phoenix/Fischar	Petroleum ether	6324035	5
fischar	Petroleum ether	6324035	1
Otto Fischar/Gehe	Petroleum ether	6324035	1
Hedinger	Petroleum ether	6324041	2
Otto Fischar/Gehe	Petroleum ether	6324044	2
Fischar	Petroleum ether	6324041	3
Caelo	Petroleum ether	6324044	1
Caelo	Petroleum ether	6324042	1
Phönix/Fischar	Petroleum ether	6324051	28
Fischar	Petroleum ether	6324055	1
Fischar/Anzag	Petroleum ether	6324063	1
Fischar	Petroleum ether	6324064	8
Fischer	Petroleum ether	6324062	1
Fischar/Sanacorp	Petroleum ether	6324065	1
Fischar	Petroleum ether	6324072	14
Fischar	Petroleum ether	6324082	15
Fischar	Petroleum ether	6324074	1
Noweda / Fischer	Petroleum ether	6324075	1
Caelo	Petroleum ether	6324082	2
Fischar/Alliance Heathc...	Petroleum ether	6324064	1
Fischar	Petroleum ether	6324075	1
Fischar	Petroleum ether	6324092	1
Fischar	Petroleum ether	6324083	1
FISCHER GMBH / GEHE	Petroleum ether	6324091	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Anzag	Petroleum ether	6324092	1
Fischar / Phoenix	Petroleum ether	6324093	1
Fischer	Petroleum ether	6324095	1
Fischer	Petroleum ether	6324102	1
Fischar	Petroleum ether	6324103	2
Fischar/Sanacorp	Petroleum ether	6324103	1
Fischer	Petroleum ether	6324103	1
Fischbar	Petroleum ether	6324103	1
Fischar	Petroleum ether	6324104	3
Fischar/Phönix	Petroleum ether	6324105	1
Fischar	Petroleum ether	6324105	4
Caelo	Petroleum ether	6324113	1
Fagron	Petroleum ether	6324114	1
Fischar	Petroleum ether	6324114	2
Otto Fischar GmbH	Petroleum ether	6324115	1
Fiebig	Petroleum ether	6324123	1
Fischar/Phönix	Petroleum ether	6324123	1
Fischar	Petroleum ether	6324123	1
Fischar	Petroleum ether	6325012	3
Fischer/Phönix	Petroleum ether	6325014	1
Fischar	Petroleum ether	6325022	2
Caelo	Petroleum ether	6325021	1
fischar	Petroleum ether	6325023	1
fischar	Petroleum ether	6325024	1
Fischar	Petroleum ether	6325023	2
Fischar	Petroleum ether	6325024	2
Fischar	Petroleum ether	6325025	1
Fischer/Phönix	Petroleum ether	6325034	3
Fischar	Petroleum ether	6325034	1
Fischar GmbH	Petroleum ether	6325044	1
Fischar	Petroleum ether	6325051	2
Fischar	Petroleum ether	6325042	1
Fischar/Phoenix	Petroleum ether	6325051	1
Fischar	Petroleum ether	6325052	3
Otto Fischar/Sanacorp	Petroleum ether	6325051	1
Fischar/Anzag	Petroleum ether	6325053	1
Caelo	Petroleum ether	6325055	2
Hubertus Apotheke	Petroleum ether	6325055	1
Otto Fischar	Petroleum ether	6325063	1
Fischar	Petroleum ether	6325063	3
Fischar/Phönix	Petroleum ether	6325081	3
Fischer	Petroleum ether	6325081	1
fischar	Petroleum ether	6325087	1
Otto Fischer / Sanacorp	Petroleum ether	6325083	1
Fischar/Phönix	Petroleum ether	6325086	1
Fischar/Gehe	Petroleum ether	6325089	1
Fischar	Petroleum ether	6325101	1
Fischar/Alliance H.	Petroleum ether	6325104	1
AEP/ Fischar	Petroleum ether	6325105	1
Fischer/Sanacorp	Petroleum ether	6325105	1
Fischer	Petroleum ether	6325105	1
Fischar	Petroleum ether	6325106	1
Fischer	Petroleum ether	6325107	1
Otto Fischer/Noweda	Petroleum ether	6325108	1
Fischar/ Sanacorp	Petroleum ether	6325108	1
Fischar/Fiebig	Petroleum ether	6325109	2
fischar	Petroleum ether	6325109	1
Frau Brinn-Holtz	Petroleum ether	6325109	1
Fischar	Petroleum ether	6325109	1
Sanacorp	Petroleum ether	632510A	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Petroleum ether	6325109	1
Caelo	Petroleum ether	6325111	2
Otto Fischar GmbH & Co...	Petroleum ether	6325111	1
Fischar Otto / Fischar ...	Petroleum ether	6325117	1
Fischar	Petroleum ether	6326011	3
Otto Fischar/Noweda	Petroleum ether	6325119	1
Fischar/Noweda	Petroleum ether	6326011	1
hedinger	Petroleum ether	6326011	1
Fischar/Phönix	Petroleum ether	6326011	1
Fischar	Petroleum ether	6326012	1
Caelo	Petroleum ether	6326014	1
fischar	Petroleum ether	6326021	2
Fischar	Petroleum ether	6326026	1
fischar/Sanacorp	Petroleum ether	6326026	1
Otto Fischer/Phönix	Petroleum ether	6326026	4
Otto Fischer/Gehe	Petroleum ether	6326027	1
Fischar/Alliance Health...	Petroleum ether	6326027	2
Alliance	Petroleum ether	6326027	1
wepa/Sanaco	Petroleum ether	6326031	1
Fischar	Petroleum ether	6326031	1
otto fischer	Petroleum ether	6326035	1
Fischar	Petroleum ether	6326032	1
Fischar/Sanacorp	Petroleum ether	6326047	1
Otto Fischer/Phönix	Petroleum ether	6326047	1
Caelo	Petroleum ether	6326055	1
Fischar/phoenix	Petroleum ether	6326061	5
Otto Fischer	Petroleum ether	6326052	1
Fischar Otto / Fischar ...	Petroleum ether	6326052	1
Fischar	Petroleum ether	6326061	1
Fischar	Petroleum ether	6326071	2
Bombastus	Petroleum ether	6326063	1
Caelo	Petroleum ether	6326063	1
Phoenix	Petroleum ether	6326072	1
Fischar	Petroleum ether	6326072	1
Fischar/Sanacorp	Petroleum ether	6326081	1
Fischar/Noweda	Petroleum ether	6326081	1
AEP/Fischar	Petroleum ether	6326092	1
Fischar	Petroleum ether	6326084	2
Fischar	Petroleum ether	6326092	1
Fischar	Petroleum ether	6326095	1
tob	Petroleum ether	6326092	1
Caelo	Petroleum ether	6326104	1
Otto Fischer	Petroleum ether	6326105	1
phönix	Petroleum ether	6326112	1
Fischar	Petroleum ether	6327013	1
Fischar	Petroleum ether	6327017	2
Fischar	Petroleum ether	6327019	1
Fischer/Phönix	Petroleum ether	6327018	1
Fischar/Sanacorp	Petroleum ether	6327018	2
Noweda	Petroleum ether	6327021	1
Fischar	Petroleum ether	6327024	1
Fischer	Petroleum ether	6327024	1
Otto Fischer/Noweda	Petroleum ether	6327024	1
Fischar	Petroleum ether	6327034	1
Fischar	Petroleum ether	6327054	1
Phönix/Fischer	Petroleum ether	6327034	1
Fischar Otto / Fischar ...	Petroleum ether	6327051	1
Fischar/ Noweda	Petroleum ether	632704	1
Fischer/Phönix	Petroleum ether	6327061	1
Gehe/Fischar	Petroleum ether	6327062	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Petroleum ether	6327062	2
Fischar	Petroleum ether	6327074	1
Phönix	Petroleum ether	6327093	1
Otto Fischar	Petroleum ether	6327072	1
phönix	Petroleum ether	6327062	1
Fischar/ AHD	Petroleum ether	6327093	1
Noweda	Petroleum ether	6327093	1
Fischar	Petroleum ether	6327102	2
Otto Fischar/Sanacorp	Petroleum ether	6327116	1
L10/Fischar	Petroleum ether	632801A	1
Sanacorp	Petroleum ether		2
Fischar	Petroleum ether	6327101	1
Noweda	Petroleum ether	6327105	1
ACM	Petroleum ether	KNC0833J15	1
Otto Fischar GmbH & Co...	Petroleum ether	X81A41	1
Hedinger/Otto	Petroleum ether		1
Hedinger/N/08.05.13/PZN...	Petroleum ether	28167B	1
Noweda	Petroleum ether	28164B	1
Gehe	Petroleum ether	28166B	1
Fischar Phönix	Petroleum ether	6323094	1
Stumpf/Hedinger	Petroleum ether	281718	1
Hedinger/Phönix	Petroleum ether	28172B	1
A.H.Hedinger vdL	Petroleum ether	28172B	1
Anzag/Hedinger	Petroleum ether	28172B	1
Hedinger/Gehe	Petroleum ether	28172b	1
Hedinger/AHD	Petroleum ether	28178B	1
Fischar Phönix	Petroleum ether	6323112	1
Hedinger/Phönix	Petroleum ether	28167B	2
Phoenix/hedinger	Petroleum ether	28169B	2
Phoenix/hedinger	Petroleum ether	28168B	1
Fischar	Petroleum ether	7283797	1
Fischar	Petroleum ether	6323033	1
Fischar	Petroleum ether	6322073	1
Hedinger/Sanacorp	Petroleum ether	28169B	4
Fischar / Gehe	Petroleum ether	6327033	1
Hedinger	Petroleum ether	204311170	1
Hedinger	Petroleum ether	(028)028	1
Hedinger/Noweda	Petroleum ether	21728	1
Hedinger / Anzag	Petroleum ether	28163B	1
Hedinger	Petroleum ether	281668	2
Hedinger/Anzag 11.06.2012	Petroleum ether	28166B	1
Anzag 21.08.2012	Petroleum ether	28166b	1
Hedinger, Noweda	Petroleum ether	28166B	2
Hedinger/Sanacorp	Petroleum ether	28166B	1
HEDINGER	Petroleum ether	28166B	2
Noweda	Petroleum ether	28166B	1
HEDINGER	Petroleum ether	28167B	1
Hedinger, Noweda	Petroleum ether	28167B	3
Caelo	Petroleum ether	28167B	1
Gehe-Hedinger	Petroleum ether	28167B	1
Kruger/Hedinger/06.01.14	Petroleum ether	28167B	1
anzag 15.11.12	Petroleum ether	28167B	1
Hedinger/Stumpf	Petroleum ether	28168B	1
Phönix	Petroleum ether	28168B	1
Hedinger, Noweda	Petroleum ether	28168B	2
Noweda/Hedinger	Petroleum ether	28168B	2
Hedinger/ Gehe	Petroleum ether	28168B	1
HEDINGER / 20	Petroleum ether	28168B	1
Hedinger / Anzag	Petroleum ether	28168B	1
Hedinger/Otto	Petroleum ether	281698	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/Sanacorp	Petroleum ether	281698	1
Hedinger/Phönix	Petroleum ether	28169B	3
Phönix	Petroleum ether	28169B	1
Stumpf / Hedinger	Petroleum ether	28169B	1
Noweda	Petroleum ether	28169B	1
Hedinger / Noweda	Petroleum ether	28169B	1
A.Hedinger / Kehr	Petroleum ether	281708	1
Dr. Seeger/ Phönix	Petroleum ether	28169B	1
Caelo	Petroleum ether	28170B	1
Noweda/Hedinger	Petroleum ether	28170B	1
Hedinger/Holdermann	Petroleum ether	28170B	1
Hedinger1	Petroleum ether	28170B	1
Hedinger	Petroleum ether	281718	1
Hedinger/V.d.L.	Petroleum ether	28171B	1
Jenne	Petroleum ether	28171B	2
Hedinger / Alliance HC	Petroleum ether	28171B	1
Hedinger/ Phönix	Petroleum ether	28171B	1
Phoenix	Petroleum ether	28171B	1
Stumpf/Hedinger	Petroleum ether	28172B	1
HEDINGER / 30	Petroleum ether	28172B	2
Hedinger / Phoenix	Petroleum ether	28172B	1
Gehe-Hedinger	Petroleum ether	28172B	1
Hedinger/ Phoenix WE. . .	Petroleum ether	28172B	1
AllianceHealthcare	Petroleum ether	28172b	1
Hedinger/Gehe	Petroleum ether	281738	1
HEDINGER	Petroleum ether	28173B	1
Hedinger/ Sanacorp	Petroleum ether	28173B	1
Hedinger/ Gehe	Petroleum ether	28173B	3
Hedinger über Sanacorp	Petroleum ether	28173B	1
Jenne	Petroleum ether	28173B	1
Sanacorp 8.12.15	Petroleum ether	28174B	2
Hedinger / Gehe	Petroleum ether	28174B	2
Caelo	Petroleum ether	281748	1
Heidinger/Anzag	Petroleum ether	28174B	1
Hedinger/Anzag	Petroleum ether	28174B	1
Noweda/Hedinger	Petroleum ether	28174B	1
Euro OTC	Petroleum ether	28174B	1
HEDINGER/Phönix	Petroleum ether	28174B	2
Hedinger/Phoenix	Petroleum ether	28174B	1
Hedinger	Petroleum ether	281758	1
Sanacorp WE:04.02.2017 ...	Petroleum ether	28175B	1
phönix	Petroleum ether	28175B	1
Hedinger/Gehe	Petroleum ether	28176B	1
Phönix EK:6,26 EUR -11	Petroleum ether	28176B	1
Hedinger/ Phoenix	Petroleum ether	28176B	1
Hedinger/Noweda	Petroleum ether	28176B	1
Hedinger	Petroleum ether	281778	1
hedinger	Petroleum ether	298171B	1
Hedinger	Petroleum ether	12121401	1
Hedinger	Petroleum ether	2B175B	1
Hedinger/Gehe	Petroleum ether	1231234	1
Caelo	Petroleum ether	13141302	1
Hedinger	Petroleum ether	13081309	1
Hedinger	Petroleum ether	15091402	1
Hedinger / Gehe	Petroleum ether	1301401	1
Hedinger, Heiligenwiese...	Petroleum ether		1
Euro OTC	Petroleum ether	2304123	1
Kost/kost	Petroleum ether	4065061213	1
Fischar/Anzag	Petroleum ether	6323016	1
Fischar	Petroleum ether	6323032	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischar/Anzag	Petroleum ether	6323024	1
Fischar/ Phönix	Petroleum ether	6323081	12
Fischar/Gehe	Petroleum ether	632309A	1
Fischar/Gehe	Petroleum ether	6324062	1
anzag	Petroleum ether	6324063	1
Anzag/Fischar	Petroleum ether	6324064	1
Fischar	Petroleum ether	6324091	1
Fischar	Petroleum ether	6324115	1
Fischar/ Noweda	Petroleum ether	6324104	1
otto fischar/Phoenix	Petroleum ether	6325044	1
Fischar/Fiebig	Petroleum ether	6325065	1
Fischar / Gehe	Petroleum ether	6325055	1
Fischar	Petroleum ether	6325065	1
Fischar	Petroleum ether	6326027	1
Fischar	Petroleum ether	6325081	1
Fischar / Gehe	Petroleum ether	6326014	1
Fischar/ Phönix	Petroleum ether	6325106	1
Fischar / Gehe	Petroleum ether	6326061	1
Fischar / Gehe	Petroleum ether	6326075	1
Fischbar	Petroleum ether	6326112	1
Caelo	Petroleum ether	6327016	1
Fischar	Petroleum ether	6326075	1
Fischar Otto / Fischar ...	Petroleum ether	6326113	1
Fischar / Gehe	Petroleum ether	6327021	1
Caelo	Petroleum ether	64/AKD53	1
Hedinger/Sanacorp	Petroleum ether		1
Großhandel	Petroleum ether		1
Alliance/ Hedinger	Petroleum ether	Y50L63	1

- 90 193 spectra from 501 *Apo-Ident* customers from a total of 26 760 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Petroleum ether* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Petroleum ether* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	140	0	35 404
Type B	0	116	4	93 754
Type C	0	1140	19	90 193

The substance/substance group *Petroleum ether* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.7143 %)
Type B	100.0000 % (> 99.9761 %)	96.6667 % (> 94.1667 %)
Type C	100.0000 % (> 99.6586 %)	98.3607 % (> 98.1018 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30899	30899	0.00	33.22
31253	31253	0.00	37.72
31283	31283	0.00	33.59

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Petroleum oil, white, cohobated
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31212-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Petroleum oil, white, cohobated; Oleum petrae album rectific.; Petroleum album rectific.

Special notes

When selecting the *Petroleum oil, white, cohobated* substance/substance group, the following information is displayed to the user:

- Advice for the identification of petroleum oil: According to the German BfArM, petroleum oil is considered as critical pharmaceutical in compliance to German § 5 AMG. Critical pharmaceuticals must not be placed on the market or be applied to humans.

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Petroleum oil, white, cohobated*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Petroleum oil, w...	13146704	31212	40	1403110
Caelo	Petroleum oil, w...	14019607	31387	40	20140124

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Petroleum oil, white, cohobated*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Petroleum oil, white, cohobated*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Petroleum oil, white, cohoba...	10356508	30411	40
Caelo	Petroleum oil, white, cohoba...	14019607	31387 [†]	20
Caelo	Petroleum oil, white, cohoba...	163278	33109	40
Caelo	Petroleum oil, white, cohoba...	182344	34238	40

- 93 734 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 22 spectra from 10 *Apo-Ident* customers from 15 batches from the substance/substance group *Petroleum oil, white, cohobated*.
- These include spectra of independent samples from 14 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Petroleum oil, white, cohoba...	13086401	1
Caelo	Petroleum oil, white, cohoba...	14242203	2
Caelo	Petroleum oil, white, cohoba...	14242201	1
Caelo	Petroleum oil, white, cohoba...	2306921	1
Caelo	Petroleum oil, white, cohoba...	11379102	1
Caelo	Petroleum oil, white, cohoba...	12173101	2
Caelo	Petroleum oil, white, cohoba...	13146702	3
Caelo	Petroleum oil, white, cohoba...	13086402	1
Caelo	Petroleum oil, white, cohoba...	13303603	2
Caelo	Petroleum oil, white, cohoba...	11198403	1
Caelo	Petroleum oil, white, cohoba...	14019607	3
Caelo	Petroleum oil, white, cohoba...	15035705	1
Caelo	Petroleum oil, white, cohoba...	14130105	1
Caelo	Petroleum oil, white, cohoba...	14019602	1
Caelo	Petroleum oil, white, cohoba...	12173103	1

- 91 330 spectra from 501 *Apo-Ident* customers from a total of 26 994 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Petroleum oil, white, cohobated* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Petroleum oil, white, cohobated* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	140	0	93 734
Type C	0	18	4	91 330

The substance/substance group *Petroleum oil, white, cohobated* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.7143 %)
Type C	100.0000 % (> 99.6590 %)	81.8182 % (> 68.1818 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31212	31212	0.00	33.79
31387	31212	4.57	33.39

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Phytomenadione**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30747-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Phytomenadione; Phytomenadionum

Special notes

When selecting the *Phytomenadione* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Phytomenadione*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Phytomenadione	12H31-N02	30747	40	1405178
Fagron	Phytomenadione	14I02-B04-301537	31621	40	20141210
Fagron	Phytomenadione	16D07-B03-328730	33304	40	20160603

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Phytomenadione*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Phytomenadione*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Phytomenadione	14I02-B04-301537	31621 [†]	20
Fagron	Phytomenadione	16D07-B03-328576	33303	40
Fagron	Phytomenadione	17F12-B12-344796	34086	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 5 spectra from 4 *Apo-Ident* customers from 5 batches from the substance/substance group *Phytomenadione*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Phytomenadione	13D26-B07-286714	1
Merck	Phytomenadione	EH14000190	1
Merck	Phytomenadione	EH13000590	1
Fagron	Phytomenadione	13C25-N07	1
Fagron	Phytomenadione	17F28-H10-00239	1

- 91 347 spectra from 501 *Apo-Ident* customers from a total of 27 004 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Phytomenadione* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Phytomenadione* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	100	0	93 774
Type C	0	5	0	91 347

The substance/substance group *Phytomenadione* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6603 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30747	30747	0.00	116.01
31621	30747	4.74	119.08
33304	30747	8.17	115.48

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Pine silvestris oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31183-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Pine silvestris oil; Oleum pini silvestris

Special notes

When selecting the *Pine silvestris oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Pine silvestris oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Pine silvestris ...	13E24-N05	31183	40	1403166
Caelo	Pine silvestris ...	13334502	31202	40	1403168
Caelo	Pine silvestris ...	13334501	31405	40	20131009

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Pine silvestris oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Pine silvestris oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Pine silvestris oil	10084705	30252	40
Caelo	Pine silvestris oil	16298401	33250	40
Caelo	Pine silvestris oil	16298401	33424	40

- 93 754 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 58 spectra from 32 *Apo-Ident* customers from 32 batches from the substance/substance group *Pine silvestris oil*.
- These include spectra of independent samples from 30 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Pine silvestris oil	13334501	4
noweda	Pine silvestris oil	13334501	2
Phönix	Pine silvestris oil	13334501	2
Caelo	Pine silvestris oil	14252202	3
Caelo	Pine silvestris oil	14252203	2
Fagron	Pine silvestris oil	14J27-B06	1
Caelo	Pine silvestris oil	15239101	4
Caelo	Pine silvestris oil	15239102	1
Fagron	Pine silvestris oil	15F12-B02-308878	1
Caelo	Pine silvestris oil	16007406	1
Taoasis	Pine silvestris oil	20017-118846	1
Taoasis	Pine silvestris oil	28663-119531	1
Fagron	Pine silvestris oil	15F12-B02-323109	1
Caelo	Pine silvestris oil	26091404	1
Fagron	Pine silvestris oil	15F12B02308878	1
Taoasis	Pine silvestris oil	28663-119029	1
Taoasis	Pine silvestris oil	28663-120876	2
Taoasis	Pine silvestris oil	35808-123038	2
Taoasis	Pine silvestris oil	49977-12442	1
Taoasis	Pine silvestris oil	49977-125240	2
Taoasis	Pine silvestris oil	8169-117244	4
Caelo	Pine silvestris oil	11024602	1
Caelo	Pine silvestris oil	11024603	1
Caelo	Pine silvestris oil	11161605	2
Noweda	Pine silvestris oil	11161605	1
Caelo	Pine silvestris oil	12080504	4
Caelo	Pine silvestris oil	12246905	1
Caelo	Pine silvestris oil	12246907	2
Fagron	Pine silvestris oil	12G20-N04	1
Caelo	Pine silvestris oil	122466904	1
Fagron	Pine silvestris oil	13E24-N05	1
Taoasis	Pine silvestris oil	30600-18490	1
Taoasis/Phönix	Pine silvestris oil	32372-110382	1
Taoasis	Pine silvestris oil	33717-112501	2
Fagron	Pine silvestris oil	12L03-N05	1

- 91 294 spectra from 501 *Apo-Ident* customers from a total of 26 977 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Pine silvestris oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Pine silvestris oil* and it was

evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	110	10	93 754
Type C	0	54	4	91 294

The substance/substance group *Pine silvestris oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	91.6667 % (> 89.1667 %)
Type C	100.0000 % (> 99.6587 %)	93.1034 % (> 87.9310 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31183	31183	0.00	21.42
31202	31202	0.00	22.02
31405	31202	2.46	21.84

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Pine tar**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31221-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Pine tar; Pix liquida

Special notes

When selecting the *Pine tar* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Pine tar*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Pine tar	13008302	31221	40	1403320
Caelo	Pine tar	16017905	33442	40	20160217
Caelo	Pine tar	180594	34085	40	20180227*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Pine tar*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 2 reference samples of the substance/substance group *Pine tar*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Pine tar	73260048	30280	40
Caelo	Pine tar	13008309	31614	60

- 93 774 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- 0 spectra from 0 *Apo-Ident* customers from 0 batches from the substance/substance group *Pine tar*.
- These include spectra of independent samples from 0 batches, of which no spectra have been used for the generation of the database.
- 91 352 spectra from 501 *Apo-Ident* customers from a total of 27 009 batches from a further 334 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Pine tar* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Pine tar* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	60	40	93 774
Type C	0	0	0	91 352

The substance/substance group *Pine tar* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	60.0000 % (> 57.0000 %)
Type C	100.0000 % (> 99.6586 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31221	31221	0.00	164.23
33442	31221	22.53	206.09
34085	34085	0.00	282.05

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Polyhexanide solution 20%**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30596-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Polyhexanide solution 20%; Polihexanidi solutio 20%; Solutio polyhexanidi 20%

Special notes

When selecting the *Polyhexanide solution 20%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Polyhexanide solution 20%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Polyhexanide sol. ...	13H05-N04	31180	40	1403256
Fagron	Polyhexanide sol. ...	13K29-B30-288935	31571	40	20131217
Fagron	Polyhexanide sol. ...	16B16-B04-322720	32733	40	20160229
Fagron	Polyhexanide sol. ...	18B16-B10-345748	34087	40	20180306

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Polyhexanide solution 20%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Polyhexanide solution 20%*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Polyhexanide solution 20%	12A04-N06	30596	40
Fagron	Polyhexanide solution 20%	13K29-B30-288935	31571 [†]	20
Fagron	Polyhexanide solution 20%	15H28-B11-313603	32276	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 216 spectra from 120 *Apo-Ident* customers from 64 batches from the substance/substance group *Polyhexanide solution 20%*.
- These include spectra of independent samples from 60 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Polyhexanide solution 20%	13351M28	2
Fargom	Polyhexanide solution 20%	116B16-B04-320854	1
Fagron	Polyhexanide solution 20%	13H05-N04	10
Fagron	Polyhexanide solution 20%	13K29-B30	1
GEHE 14.04.2014	Polyhexanide solution 20%	13K29-B30288934	1
Fagron	Polyhexanide solution 20%	13K29-B30288934	1
Fagron	Polyhexanide solution 20%	13K29-B30-288934	5
Fagron	Polyhexanide solution 20%	13K29-B30-288938	6
Caelo	Polyhexanide solution 20%	13K29-B30-288934	1
Caelo	Polyhexanide solution 20%	13K29-B30-288938	2
Fagron	Polyhexanide solution 20%	14F25-B06	11
Sanacorp	Polyhexanide solution 20%	14f25-b06-2962291	1
Caelo	Polyhexanide solution 20%	14F25-B06-296291	2
Fagron	Polyhexanide solution 20%	14F25-B06-296291	2
Caelo	Polyhexanide solution 20%	14F25-B06-296292	2
Fagron	Polyhexanide solution 20%	14F25-B06-296292	11
Gehe 16.12.14	Polyhexanide solution 20%	14F25-B06-296292	1
Fagron	Polyhexanide solution 20%	14F25-b06-296292	1
Fragon/Phönix	Polyhexanide solution 20%	14F25-B06-300092	1
Fagron	Polyhexanide solution 20%	14F25-B06-300092	6
Fagron	Polyhexanide solution 20%	14F25-B06-300093	7
Caelo	Polyhexanide solution 20%	14F25-B06-300093	2
Caelo	Polyhexanide solution 20%	14F25-B06-305013	1
Fagron	Polyhexanide solution 20%	14F25-B06-305013	3
Noweda	Polyhexanide solution 20%	14F25B06305013	1
Fagron	Polyhexanide solution 20%	14GR132285	2
Caelo	Polyhexanide solution 20%	14GR132285	1
Fagron	Polyhexanide solution 20%	15C20-B06-309678	4
Fagron	Polyhexanide solution 20%	14GR274618	1
Fagron	Polyhexanide solution 20%	15H28-B11	6
Noweda	Polyhexanide solution 20%	15H28-B11-312309	2
Fargon/ Gehe	Polyhexanide solution 20%	15H28-B11-312309	1
Fagron	Polyhexanide solution 20%	15H28-B11-312309	2
Fagron	Polyhexanide solution 20%	15h28-b11-312309	1
Fagron	Polyhexanide solution 20%	15H28-B11-313603	7
Fagron	Polyhexanide solution 20%	15H28-B11-313604	3
Fagron	Polyhexanide solution 20%	15H28-B11313603	2
Euro OTC	Polyhexanide solution 20%	15H28-B11-313604	1
Fagron	Polyhexanide solution 20%	15H28-B11-315703	4
Fagron	Polyhexanide solution 20%	16B16-B04	6
Fagron	Polyhexanide solution 20%	16B16-b04	1
Phoenix	Polyhexanide solution 20%	16B16-B04	1
Fagron	Polyhexanide solution 20%	16B16-B04-318789	4
Caelo	Polyhexanide solution 20%	16B16-B04-318789	1
Fagron	Polyhexanide solution 20%	16B16-B04-320854	6
Fagron	Polyhexanide solution 20%	16B16-B04-322720	4
Fagron	Polyhexanide solution 20%	16B16-B04-325487	7
Fagron	Polyhexanide solution 20%	16G29-B01	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Polyhexanide solution 20%	16G259-B01-328171	1
Fagron	Polyhexanide solution 20%	16G29-B01-328171	3
Fragon/Noweda	Polyhexanide solution 20%	16G29-B01-328171	1
Ebert-Jacobi	Polyhexanide solution 20%	16G29-B01-328171	1
Fagron	Polyhexanide solution 20%	16G29-B01-330629	1
Ebert-Jacobi	Polyhexanide solution 20%	16G29-B01-330629	1
Fagron	Polyhexanide solution 20%	16G29-B01-330630	3
Fagron	Polyhexanide solution 20%	16G29-B01-335336	4
Fagron	Polyhexanide solution 20%	16g29-b01-335336	1
Fagron	Polyhexanide solution 20%	16G29-B01-337612	3
Fragon GmbH	Polyhexanide solution 20%	16GR025056	1
Fagron	Polyhexanide solution 20%	17L14-B02-343991	1
Fagron	Polyhexanide solution 20%	18071501	1
Beiersdorf	Polyhexanide solution 20%	18816-810-345748	1
Fagron	Polyhexanide solution 20%	18B10-345748	1
Fagron	Polyhexanide solution 20%	18B16-B10-345748	1
Fagron	Polyhexanide solution 20%	2215I-03580	1
Fagron	Polyhexanide solution 20%	3714M-03580	1
Fagron	Polyhexanide solution 20%	11GR180545	2
Caelo	Polyhexanide solution 20%	AKD17	1
Fagron	Polyhexanide solution 20%	12A04-N06	2
Fagron	Polyhexanide solution 20%	12K20-N01	5
Fagron	Polyhexanide solution 20%	12k20-n01	2
Fragon	Polyhexanide solution 20%	12K20-N01	1
Caelo	Polyhexanide solution 20%	12GR227646	1
Caelo	Polyhexanide solution 20%	12B09-N01	1
Noweda	Polyhexanide solution 20%	12K20-N01	1
Gehe 16.04.2013	Polyhexanide solution 20%	12K20-N01	1
Caelo	Polyhexanide solution 20%	12K20N01	1
Phoenix	Polyhexanide solution 20%	12K20-N01	1
Fagron	Polyhexanide solution 20%	13H05-N01	8
B.Braun/ Noweda	Polyhexanide solution 20%	13103M18	3
Fagron	Polyhexanide solution 20%	13H05-N040339/904	1
Fagron	Polyhexanide solution 20%	13H05-No4	1
Phytocomm	Polyhexanide solution 20%	27803278	1
Caelo	Polyhexanide solution 20%	15H28-B11-315703	1
Fagron	Polyhexanide solution 20%	12D27-N07	1
Caelo	Polyhexanide solution 20%	14F25-B06	1
Ichthyol Gesellschaft/AHD	Polyhexanide solution 20%	14F25-B06-296292	1

- 91 136 spectra from 501 *Apo-Ident* customers from a total of 26 947 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Polyhexanide solution 20%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Polyhexanide solution 20%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	100	0	93 774
Type C	0	190	26	91 136

The substance/substance group *Polyhexanide solution 20%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6586 %)	87.9630 % (> 86.5741 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31180	31180	0.00	29.33
31571	31180	8.10	31.33
32733	31180	12.18	32.99
34087	31180	4.97	26.57

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Polyol fatty acid ester
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30755-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Polyol fatty acid ester; Cetiol® HE

Special notes

When selecting the *Polyol fatty acid ester* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Polyol fatty acid ester*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Polyol fatty aci...	11162809	30755	40	20110616
Caelo	Polyol fatty aci...	13292302	31179	40	1402589
Caelo	Polyol fatty aci...	13357702	31366	40	20131112
Caelo	Polyol fatty aci...	16254506	33460	40	20160901

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Polyol fatty acid ester*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 160 spectra from 4 reference samples of the substance/substance group *Polyol fatty acid ester*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Polyol fatty acid ester	zu erfassen	30567	40
Caelo	Polyol fatty acid ester	161120	32196	40
Caelo	Polyol fatty acid ester	16254506	33459	40
Caelo	Polyol fatty acid ester	172830	33676	40

- 93 714 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 65 spectra from 43 *Apo-Ident* customers from 25 batches from the substance/substance group *Polyol fatty acid ester*.
- These include spectra of independent samples from 20 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Sanacorp	Polyol fatty acid ester	1/AKD14	1
Caelo	Polyol fatty acid ester	13357703	4
Phönix	Polyol fatty acid ester	14033502	1
Caelo	Polyol fatty acid ester	14033502	7
Caelo	Polyol fatty acid ester	143033502	1
Caelo	Polyol fatty acid ester	14322405	3
Caelo	Polyol fatty acid ester	15298004	2
Caelo	Polyol fatty acid ester	14322407	1
Caelo	Polyol fatty acid ester	16254504	2
Caelo	Polyol fatty acid ester	16254506	1
Caelo	Polyol fatty acid ester	16254502	1
Caelo	Polyol fatty acid ester	16254508	1
Caelo	Polyol fatty acid ester		1
Caelo	Polyol fatty acid ester	17072203	1
Caelo	Polyol fatty acid ester	11081403	1
Caelo	Polyol fatty acid ester	11162809	1
GEHE	Polyol fatty acid ester	12361804	2
Caelo	Polyol fatty acid ester	11162811	1
Caelo	Polyol fatty acid ester	12361804	4
Caelo	Polyol fatty acid ester	12361808	4
Caelo	Polyol fatty acid ester	13158901	1
Caelo	Polyol fatty acid ester	13158902	5
Caelo	Polyol fatty acid ester	13158903	5
Noweda	Polyol fatty acid ester	13158903	1
Caelo	Polyol fatty acid ester	13292302	4
Caelo	Polyol fatty acid ester	13158904	1
Phönix	Polyol fatty acid ester	13357702	1
Caelo	Polyol fatty acid ester	13357702	6
Caesar/ Noweda	Polyol fatty acid ester	2613E-01243	1

- 91 287 spectra from 501 *Apo-Ident* customers from a total of 26 984 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Polyol fatty acid ester* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Polyol fatty acid ester* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	160	0	93 714
Type C	0	65	0	91 287

The substance/substance group *Polyol fatty acid ester* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.2500 %)
Type C	100.0000 % (> 99.6587 %)	100.0000 % (> 90.7692 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30755	31179	8.04	17.45
31179	31179	0.00	18.93
31366	31179	5.89	19.02
33460	31179	7.39	16.65

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Polysorbate 80 / Macrogolglycerolricinoleat
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30766-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Polysorbate 80 / Macrogolglycerolricinoleat; Macrogolglycerol ricinoleate; Macrogolglyceroli ricinoleas; Poloxyaethylenglycerinum triricinoleinicum 35; Polysorbate 80; Polysorbatum 80

Special notes

When selecting the *Polysorbate 80 / Macrogolglycerolricinoleat* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Polysorbate 80 / Macrogolglycerolricinoleat*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Macrogolglycerol...	12E29-N02	30766	60	20120612
Caelo	Polysorbate 80	12116009	30893	80	20120427
Fagron	Macrogolglycerol...	12H08-N01	31015	40	20120827
Caelo	Polysorbate 80	13372015	31325	40	AR-15-FG-008176-01
Caelo	Macrogolglycerol...	13272708	31368	40	1505503
Fagron	Polysorbate 80	17A23-B09-331594	33300	40	20170208

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 300 spectra of 6 reference samples from the substance/substance group *Polysorbate 80 / Macrogolglycerolricinoleat*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 6 different batches.
- 35 244 spectra from a total of 834 batches from further 210 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 250 spectra from 7 reference samples of the substance/substance group *Polysorbate 80 / Macrogolglycerolricinoleat*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Polysorbate 80	11179606	30395	40
Caelo	Polysorbate 80	12116009	30893 [†]	10
Caelo	Macrogolglycerol ricinoleate	161303	32297	40
Caelo	Macrogolglycerol ricinoleate	161545	32338	40
Caelo	Polysorbate 80	16260201	33296	40
Caelo	Macrogolglycerol ricinoleate	172238	33707	40
Caelo	Polysorbate 80	173048	33887	40

- 93 624 spectra from a total of 2160 batches from further 435 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 300 spectra from 148 *Apo-Ident* customers from 94 batches from the substance/substance group *Polysorbate 80* / *Macrogolglycerolricinoleat*.
- These include spectra of independent samples from 90 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Macrogolglycerol ricinoleate	14195501	1
Fagron	Macrogolglycerol ricinoleate	14L08-B04-302328	2
Caelo	Macrogolglycerol ricinoleate	15312304	1
Caelo	Macrogolglycerol ricinoleate	15312307	2
Caelo	Macrogolglycerol ricinoleate	16130304	2
Caelo	Macrogolglycerol ricinoleate	16154501	1
Caelo	Macrogolglycerol ricinoleate	16154504	1
Fagron	Macrogolglycerol ricinoleate	13B08-N03	2
BASF Chem trade GmbH/Ph...	Macrogolglycerol ricinoleate	13B08-N03	1
Fagron	Macrogolglycerol ricinoleate	44242575L0	1
Fagron	Macrogolglycerol ricinoleate	13E13-N11	2
Fagron	Macrogolglycerol ricinoleate	13J24-B03-287433	1
Fagron	Macrogolglycerol ricinoleate	13J24-B03-296052	1
Fagron	Macrogolglycerol ricinoleate	12G02-N11	1
Caelo	Macrogolglycerol ricinoleate	44242575L0	1
Fagron	Macrogolglycerol ricinoleate	13B08-No3	1
Fagron	Macrogolglycerol ricinoleate	16B08-B01-322893	1
Phönix	Macrogolglycerol ricinoleate	16B08-B01-328228	1
Fagron	Macrogolglycerol ricinoleate	17J10-B02-340819	1
Phoenix	Macrogolglycerol ricinoleate	73353136W0	1
Caelo	Polysorbate 80	515I-02527	1
Fagron	Polysorbate 80	813I-02527	1
Caelo	Polysorbate 80	12315324	1
Caelo	Polysorbate 80	12315317	10
VDL 17.04.14	Polysorbate 80		1
Caelo	Polysorbate 80	12315328	14
Caelo	Polysorbate 80	13152216	13
Caelo	Polysorbate 80	13152218	1
Caelo	Polysorbate 80	13152223	11
Nachprüfung 3008	Polysorbate 80	13152223	1
Caelo	Polysorbate 80	13372012	11
Caelo	Polysorbate 80	13372015	1
Gehe 01.07.2014	Polysorbate 80	13372015	1
Caelo	Polysorbate 80	1419026	1
Caelo	Polysorbate 80	14190306	5
ICHTHYOL/ahp	Polysorbate 80	14190310	1
Caelo	Polysorbate 80	14190316	19
Gehe	Polysorbate 80	14190316	1
Caelo	Polysorbate 80	14190317	1
Fagron	Polysorbate 80	14C24-B07	1
Fagron	Polysorbate 80	14C24-B07-296371	3
Fagron	Polysorbate 80	14K27-B01-307106	2
Fargon/Kehr	Polysorbate 80	14K27-B01-307106	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Devesa/ Phönix	Polysorbate 80	14K27-Bo1-307106	1
Caelo	Polysorbate 80	15104103	10
VDL Nachliefer; 2.7.15;...	Polysorbate 80	15104103	1
Caelo	Polysorbate 80	15104112	12
Ebert u. Jacobi 24.06...	Polysorbate 80	15104112	1
Caelo	Polysorbate 80	15104113	1
Caelo	Polysorbate 80	15104121	9
Caelo	Polysorbate 80	15104126	6
Caelo	Polysorbate 80	15121404	1
Caelo	Polysorbate 80	16053402	1
Caelo	Polysorbate 80	16059402	23
Alte Apotheke	Polysorbate 80	16059402	1
Cealo/AHD	Polysorbate 80	16059402	1
Sanacorp, 4,52, 18.1.17	Polysorbate 80	16059402	1
Caelo	Polysorbate 80	16059403	1
Caelo	Polysorbate 80	1612B08327237	1
Caelo	Polysorbate 80	16260201	9
Caesar & Loretz/Noweda	Polysorbate 80	16260201	1
Caelo	Polysorbate 80	16260202	1
Caelo	Polysorbate 80	16260210	4
Cealo	Polysorbate 80	16260210	1
Fagron	Polysorbate 80	16I12-B02-327237	1
GEHE	Polysorbate 80	16B17-B09	1
Caelo	Polysorbate 80	17064705	1
Caelo	Polysorbate 80	18081502	1
Caelo	Polysorbate 80	42175394	2
Fagron	Polysorbate 80	226953-0002	1
Caelo	Polysorbate 80	25062015C	1
Fagron	Polysorbate 80	808ND3713	1
Fagron	Polysorbate 80	908ND3714	2
Caelo	Polysorbate 80	11179606	1
Caelo	Polysorbate 80	10235202	1
Fargon	Polysorbate 80	10114N07	1
Caelo	Polysorbate 80	11179618	1
Bombastus	Polysorbate 80	11347907	1
Caelo	Polysorbate 80	11347907	3
Caelo	Polysorbate 80	11347909	1
Phoenix	Polysorbate 80	12116005	1
Fagron	Polysorbate 80	11L22-N12	1
Caelo	Polysorbate 80	12116005	4
Caelo	Polysorbate 80	12116009	8
Caelo	Polysorbate 80	12116013	7
Phönix 04.07.2013	Polysorbate 80	12116013	1
Caelo	Polysorbate 80	12315308	10
Caelo	Polysorbate 80	123125308	1
Caelo	Polysorbate 80	12315305	1
Caelo	Polysorbate 80	12116014	1
Alliance Healthcare	Polysorbate 80	12315308	1
phönix	Polysorbate 80	12315328	1
Fagron	Polysorbate 80	12K05-N04	4
gehe 28.12.2013	Polysorbate 80	13152216	3
Caelo	Polysorbate 80	13152219	1
Krieger	Polysorbate 80	131522232x	1
Noweda	Polysorbate 80	13A10-N03	1
Fagron	Polysorbate 80	13A10-N03	2
Fagron	Polysorbate 80	13C19-N01	1
Fagron	Polysorbate 80	13C19N01	1
Caelo	Polysorbate 80	20210312-2	1
Caelo	Polysorbate 80	21061301	1
Caelo	Polysorbate 80	27111308	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Polysorbate 80	7/AKD20	1
Caelo	Polysorbate 80	2772013	1
Caelo	Polysorbate 80	Ch12116005	1
Fagron	Polysorbate 80	12D02-N02	1
Gehe	Polysorbate 80	12K05-N04	1
Caelo	Polysorbate 80	2915A-02527	1
Klenk	Polysorbate 80	12315328	1
Caelo	Polysorbate 80	3212I-02527	1
Caelo	Polysorbate 80	3213E2527	1
Caelo	Polysorbate 80	4512A-02527	1
Caelo	Polysorbate 80	4513M	1
Caelo	Polysorbate 80		1

- 91 052 spectra from 501 *Apo-Ident* customers from a total of 26 917 batches from a further 332 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Polysorbate 80 / Macrogolglycerolricinoleat* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Polysorbate 80 / Macrogolglycerolricinoleat* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	300	0	35 244
Type B	0	250	0	93 624
Type C	0	299	1	91 052

The substance/substance group *Polysorbate 80 / Macrogolglycerolricinoleat* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.6000 %)
Type C	100.0000 % (> 99.6586 %)	99.6667 % (> 98.6667 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30766	31368	3.37	23.76
30893	31325	4.98	13.11
31015	31368	5.58	25.82
31325	31325	0.00	12.75
31368	31368	0.00	23.92
33300	31325	5.57	13.82

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Potash soap DAC**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30987-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Potash soap DAC; Sapo kalinus DAC

Special notes

When selecting the *Potash soap DAC* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Potash soap DAC*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Potash soap DAC	13008101	30987	40	1402565
Caelo	Potash soap DAC	14143401	31480	40	20140520
Caelo	Potash soap DAC	14240801	31561	40	20140815
Caelo	Potash soap DAC	14289101	31650	40	20141017

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Potash soap DAC*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 440 spectra from 12 reference samples of the substance/substance group *Potash soap DAC*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Potash soap DAC	11136402	30336	40
Caelo	Potash soap DAC	13322201	31125	40
Caelo	Potash soap DAC	14240801	31561 [†]	20
Caelo	Potash soap DAC	14289101	31650 [†]	20
Caelo	Potash soap DAC	160908	32392	40
Caelo	Potash soap DAC	160909	32393	40
Caelo	Potash soap DAC	160910	32394	40
Caelo	Potash soap DAC	171575	33408	40
Caelo	Potash soap DAC	170461	33409	40
Caelo	Potash soap DAC	170462	33410	40
Caelo	Potash soap DAC	173023	34016	40
Caelo	Potash soap DAC	180707	34149	40

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 434 spectra from a total of 2155 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 39 spectra from 20 *Apo-Ident* customers from 29 batches from the substance/substance group *Potash soap DAC*.
- These include spectra of independent samples from 27 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Potash soap DAC	11198901	1
Caelo	Potash soap DAC	12060401	1
Caelo	Potash soap DAC	11276401	1
Caelo	Potash soap DAC	12060701	1
Caelo	Potash soap DAC	13161401	2
Caelo	Potash soap DAC	13161901	2
Caelo	Potash soap DAC	13322201	1
Caelo	Potash soap DAC	13322601	1
Caelo	Potash soap DAC	13322401	1
Caelo	Potash soap DAC	14032801	1
Caelo	Potash soap DAC	14032901	1
Caelo	Potash soap DAC	14033001	1
Klenk	Potash soap DAC	14143401	1
Caelo	Potash soap DAC	14240701	1
Caelo	Potash soap DAC	14143601	1
Caelo	Potash soap DAC	14240901	1
Caelo	Potash soap DAC	14241001	2
Caelo	Potash soap DAC	14241101	2
Caelo	Potash soap DAC	14289101	2
Caelo	Potash soap DAC	14289201	1
Caelo	Potash soap DAC	15240501	4
Caelo	Potash soap DAC	15/134/KS	1
Caelo	Potash soap DAC	15240401	1
Caelo	Potash soap DAC	15435201	2
Caelo	Potash soap DAC	15182001	1
Caelo	Potash soap DAC	16091101	1
Caelo	Potash soap DAC	16091201	1
Caelo	Potash soap DAC	16256901	1
Caelo	Potash soap DAC	17221701	2

- 91 313 spectra from 501 *Apo-Ident* customers from a total of 26 980 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Potash soap DAC* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Potash soap DAC* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	371	69	93 434
Type C	0	35	4	91 313

The substance/substance group *Potash soap DAC* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	84.3182 % (> 83.6364 %)
Type C	100.0000 % (> 99.6588 %)	89.7436 % (> 82.0513 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30987	30987	0.00	58.87
31480	30987	11.87	56.51
31561	30987	12.90	58.07
31650	30987	16.34	59.98

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Propylene glycol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30872-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Propylene glycol; Propylenglycolum

Special notes

When selecting the *Propylene glycol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Propylene glycol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Propylene glycol	12372405	30872	40	20130107
Henry Lamotte	Propylene glycol	0123095	30920	40	1502554
Caelo	Propylene glycol	14346604	31663	40	1508143
Euro OTC	Propylene glycol	1509031	31738	45	20151005*
Caelo	Propylene glycol	15335403	32074	40	20150930
Caelo	Propylene glycol	15355108	32486	40	20151021
Bombastus	Propylene glycol	301162	32728	40	20160323

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 285 spectra of 7 reference samples from the substance/substance group *Propylene glycol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 7 different batches.
- 35 259 spectra from a total of 833 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 315 spectra from 9 reference samples of the substance/substance group *Propylene glycol*.
- These include spectra of independent samples from 7 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Propylene glycol	14346604	31663 [†]	20
Euro OTC	Propylene glycol	1509031	31738 [†]	15
Euro OTC	Propylene glycol	1603034	32190	40
Euro OTC	Propylene glycol	1606014	32466	40
Euro OTC	Propylene glycol	1610031	32836	40
Euro OTC	Propylene glycol	1703024	33447	40
Euro OTC	Propylene glycol	1708002	33665	40
Euro OTC	Propylene glycol	1711023	33923	40
Euro OTC	Propylene glycol	1803023	34093	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93559 spectra from a total of 2158 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2743 spectra from 428 *Apo-Ident* customers from 454 batches from the substance/substance group *Propylene glycol*.
- These include spectra of independent samples from 443 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Propylene glycol	16C10-B04-319299	2
Fagron	Propylene glycol	12J04-N06	13
Caelo	Propylene glycol	12169803	12
Sanacorp	Propylene glycol	123095	3
Caelo	Propylene glycol	1071502	1
Caelo	Propylene glycol	1140817	1
Henry Lamotte Oils GmbH...	Propylene glycol	123095	1
Henry Lamotte/Phoenix	Propylene glycol	123095	11
Großhandel	Propylene glycol	123095	23
HenryLamotteOilsGmbH	Propylene glycol	123095	1
Caelo	Propylene glycol	123095	1
Lamotte/ Finze	Propylene glycol	123095	8
Lamotte Oils GmbH	Propylene glycol	123095	1
Caelo	Propylene glycol	125099	16
Henry Lamotte Oils GmbH...	Propylene glycol	123095	1
Henry Lamotte oils GmbH	Propylene glycol	125099	1
Alliance Healthcare/Hen...	Propylene glycol	125099	2
Euro OTC	Propylene glycol	125099	9
Sanacorp	Propylene glycol	125099	5
Henry Lamotte / Allianc...	Propylene glycol	125099	1
Henry Lamotte	Propylene glycol	125099	2
Henry/Geibo	Propylene glycol	125099	1
Henry Lamotte Oils / No...	Propylene glycol	125099	2
Henry Lamotte Oils GmbH	Propylene glycol	125099	4
Henry Lamotte Oils GmbH...	Propylene glycol	125099	1
Henry Lamotte/Noweda	Propylene glycol	125099	2
Henry Lamotte GmbH	Propylene glycol	125099	1
Lamotte GmbH	Propylene glycol	125099	1
Fagron	Propylene glycol	2021621	1
Caelo	Propylene glycol	2310070	1
Phoenix	Propylene glycol	125099	1
Hedinger	Propylene glycol	29702B	1
Hedinger/ Sanacorp	Propylene glycol	29743B	1
Hedinger/Phoenix	Propylene glycol	29762B	2
Phoenix	Propylene glycol	29722B	1
Hedinger	Propylene glycol	29797B	4
Hedinger/Alliance Heath...	Propylene glycol	29775B	1
Hedinger	Propylene glycol	29775B	1
Caelo	Propylene glycol	29814B	1
Hedinger/Sanacorp	Propylene glycol	29823B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger	Propylene glycol	29837B	1
Hedinger/VDL	Propylene glycol	29847B	1
Hedinger	Propylene glycol	29823B	1
Hedinger	Propylene glycol	29859B	1
Hedinger/VDL	Propylene glycol	29859B	1
Hedinger	Propylene glycol	29859BP	1
Caelo	Propylene glycol	29879B	1
Hedinger	Propylene glycol	29905B	1
Caelo	Propylene glycol	3071402	1
Henry Lamotte	Propylene glycol	2992	1
Hedinger	Propylene glycol	29897B	1
Fagron	Propylene glycol	5021520	2
Fagron	Propylene glycol	50515PJK01	1
Caelo	Propylene glycol	514E-02110	1
Fagron	Propylene glycol	6071520	3
Fagron	Propylene glycol	8021621	2
Caelo	Propylene glycol	6071520	1
Henry Lamotte Olis GmbH...	Propylene glycol	10000001132	1
Caelo	Propylene glycol	8051503	1
Henry Lamotte	Propylene glycol	1000001132	13
Lamotte/Anzag11	Propylene glycol	1000001132	3
Caelo	Propylene glycol	1000001132	5
Henry Lamotte / Allianc...	Propylene glycol	1000001132	2
Henry Lamotte Oils GmbH	Propylene glycol	1000001132	4
Lamotte / AHD	Propylene glycol	1000001132	1
Henry Lamotte Oils / No...	Propylene glycol	1000001132	1
Henry Lamotte/ Noweda	Propylene glycol	1000001132	2
Phoenix	Propylene glycol	1000001132	1
Lamotte Sanacorp	Propylene glycol	1000001132	1
Lamotte	Propylene glycol	1000001132	1
Caelo	Propylene glycol	1000002768	1
Lamotte / Noweda	Propylene glycol	1000002768	1
Lamotte	Propylene glycol	1000002768	2
Henry Lamotte / Noweda	Propylene glycol	1000002768	1
Henr	Propylene glycol	1000002768	1
Henry Lamotte	Propylene glycol	1000002768	3
Caelo	Propylene glycol	1000002992	2
Henry Lamotte	Propylene glycol	1000002992	4
Lamotte	Propylene glycol	1000002992	1
Caelo	Propylene glycol	1113102	1
Henry Lamotte/Sanacorp	Propylene glycol	100001132	1
Phönix	Propylene glycol	114113102	1
Fagron	Propylene glycol	1304C009	6
Caelo	Propylene glycol	12071403	1
Fagron	Propylene glycol	1304C009-281100	3
Noweda 29.09.2014	Propylene glycol	1304C010-281099	1
Caelo	Propylene glycol	12345	1
Euro OTC	Propylene glycol	1307021-01	10
Fagron	Propylene glycol	1307021-01	1
Caelo	Propylene glycol	13256704	29
Caelo	Propylene glycol	13111405	1
Fagron	Propylene glycol	1314Q2110	1
Caelo	Propylene glycol	13256705	18
Noweda	Propylene glycol	13256705	3
Caesar & Loretz GmbH Hi...	Propylene glycol	13256705	1
GEHE 06.05.2014	Propylene glycol	13256705	2
Caelo	Propylene glycol	13256706	22
Sanacorp	Propylene glycol	13256706	1
Caelo	Propylene glycol	13256707	29
Caelo	Propylene glycol	13256708	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Cealo / Noweda	Propylene glycol	13256708	2
Cealo	Propylene glycol	13256709	14
Fagron	Propylene glycol	13E22-B03	11
VDL	Propylene glycol	13256709	1
Phoenix am 10.03.2014	Propylene glycol	13256709	1
Cealo	Propylene glycol	13E22-B03	1
Fagron	Propylene glycol	13E22-B03-283450	24
Phoenix	Propylene glycol	13E22B03283450	1
Fagron	Propylene glycol	13E22-B03-28350	1
Fagron	Propylene glycol	13G09-B01	5
Fagron	Propylene glycol	13G09B01	2
Fagron	Propylene glycol	13G09-B01-283456	23
Cealo	Propylene glycol	13G09-B01283456	1
Fagron	Propylene glycol	13G09B01283456	1
EuRho/Noweda	Propylene glycol	140020-01	1
Euro OTC	Propylene glycol	1401010-01	8
Cealo	Propylene glycol	140101001	2
gehe	Propylene glycol	1401010-01	2
Euro OTC	Propylene glycol	140500301	3
Fagron	Propylene glycol	1405003-01	1
Cealo	Propylene glycol	1405003-01	1
Euro OTC	Propylene glycol	1405003-01	2
Phönix 4.3.15	Propylene glycol	1405003-01	1
Cealo	Propylene glycol	14051802	21
Cealo	Propylene glycol	140510804	1
Gehe	Propylene glycol	14051802	1
noweda	Propylene glycol	14051802	1
Cealo	Propylene glycol	14051803	8
Caesar/Noweda	Propylene glycol	14051803	2
Cealo	Propylene glycol	14051804	14
VDL 14.05.14; Preis 3.07	Propylene glycol		1
VDL 14.05.14; Preis 3.07	Propylene glycol	14051804	2
Caesar & Loretz GmbH	Propylene glycol	14051804	2
Cealo	Propylene glycol	14051805	23
Phoenix	Propylene glycol	14051805	1
Euro OTC	Propylene glycol	1408020-01	12
BASF AG	Propylene glycol	140802001	1
Euro OTC	Propylene glycol	140802001	1
Cealo	Propylene glycol	14113102	18
Phönix	Propylene glycol	141113102	1
Phönix	Propylene glycol	14113102	1
Sanacorp	Propylene glycol	14113103	1
Cealo	Propylene glycol	14113103	5
Cealo	Propylene glycol	14113105	12
VDL	Propylene glycol	14113103	1
Alliance 05.2.15	Propylene glycol	14113105	1
Cealo	Propylene glycol	1424870	1
Cealo	Propylene glycol	14248702	12
Cealo	Propylene glycol	14248703	19
Sanacorp	Propylene glycol	14248703	1
Cealo	Propylene glycol	14248704	7
Euro/ Sanacorp	Propylene glycol	14248704	1
Caesar & Loretz GmbH Hi...	Propylene glycol	14248706	1
Cealo	Propylene glycol	14248706	5
Cealo	Propylene glycol	14248707	12
Sanacorp	Propylene glycol	14248706	1
Caesar & Loretz GmbH Hi...	Propylene glycol	14248707	1
AHD 26.01.15	Propylene glycol	14248707	1
Cealo	Propylene glycol	14346602	4
Euro OTC	Propylene glycol	14346602	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Propylene glycol	1428706	1
Caelo	Propylene glycol	1428704	2
caesar	Propylene glycol	14248707	1
Gehe	Propylene glycol	14346602	1
Caelo	Propylene glycol	14346603	11
Caelo	Propylene glycol	14346604	9
Caelo	Propylene glycol	14346605	2
Caelo	Propylene glycol	14346606	9
Caelo	Propylene glycol	14346607	21
Sanacorp	Propylene glycol	14346607	1
Caelo	Propylene glycol	14346608	12
Caelo	Propylene glycol	14346609	17
Caelo	Propylene glycol	1436604	1
Fagron	Propylene glycol	14A09-B04	15
Fagron	Propylene glycol	14A09B04	4
Fagron	Propylene glycol	14A09-B04-289597	21
Gehe	Propylene glycol	14A09-B04-289597	2
Phoenix	Propylene glycol	14a09-b04-289597	1
Phoenix	Propylene glycol	14A09-B04-289597	1
Caelo	Propylene glycol	14A09-B04-289598	1
Fagron	Propylene glycol	14A09-B04-289598	4
Fagron	Propylene glycol	14A09-B04-289599	4
Fagron	Propylene glycol	14C05-B02	4
Fagron	Propylene glycol	14A09-B04-289600	1
Fagron	Propylene glycol	14C05-B02-295787	10
Phönix Gotha	Propylene glycol	14C05-B02-295787	1
Gehe	Propylene glycol	14C05-B02-295787	1
Fargron	Propylene glycol	14C05-B02-295787	2
Fagron	Propylene glycol	14C05B02295787	2
Gehe 3.3.15	Propylene glycol	14C05-B02-295787	1
Fagron	Propylene glycol	14C05-B02-295788	4
Fagron	Propylene glycol	14C09-B02-289597	1
Fagron	Propylene glycol	14H22-B02	1
Fagron	Propylene glycol	14H22-B02-302506	5
BASF	Propylene glycol	1501013-01	2
BASF AG.	Propylene glycol	1501013-01	1
Euro OTC	Propylene glycol	1501013-01	15
Euro OTC	Propylene glycol	150101301	2
Linden-Apo	Propylene glycol	1501013-01	1
Euro OTC	Propylene glycol	1504031-01	21
Euro OTC	Propylene glycol	150403101	3
Euro OTC	Propylene glycol	15040301-01	1
EuRho/Noweda	Propylene glycol	1504031-01	1
Euro PTC Pharma GmbH/Sa...	Propylene glycol	1504031-01	1
Caelo	Propylene glycol	15069201	5
Sanacorp	Propylene glycol	15069201	1
Caelo	Propylene glycol	15069202	17
Noweda	Propylene glycol	15069202	1
Caelo	Propylene glycol	15081501	1
Euro OTC	Propylene glycol	1509031-01	19
EuRho	Propylene glycol	1509031-01	2
Caelo	Propylene glycol	150903101	1
Euro OTC	Propylene glycol	150903101	1
Euro OTC	Propylene glycol	1509031-02	1
Caelo	Propylene glycol	15181601	20
Euro OTC	Propylene glycol	1509031-031	1
Caelo	Propylene glycol	15181602	12
VDL Nachliefer; 11.7.15...	Propylene glycol	15181602	1
Phoenix	Propylene glycol	15181602	1
Caelo	Propylene glycol	15181603	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Propylene glycol	15181604	8
Caesar&Loretz GmbH/Anzag	Propylene glycol	15181605	2
Caelo	Propylene glycol	15181605	11
Caelo	Propylene glycol	151816055	1
Caelo	Propylene glycol	1518601	1
Caelo	Propylene glycol	1518602	1
Caelo	Propylene glycol	15262401	16
Caelo	Propylene glycol	15262402	9
Noweda	Propylene glycol	15262402	1
Caelo	Propylene glycol	15262403	21
Euro OTC	Propylene glycol	15262403	1
Caesar & Loretz GmbH	Propylene glycol	15262403	1
Caelo	Propylene glycol	15262404	14
Caelo	Propylene glycol	15262406	5
Caelo	Propylene glycol	15262407	9
Fagron	Propylene glycol	15262407	1
Caelo	Propylene glycol	15262408	9
Caelo	Propylene glycol	15262409	10
Noweda	Propylene glycol	15262409	1
Caelo	Propylene glycol	15335405	1
Caelo	Propylene glycol	1526409	1
Caelo	Propylene glycol	1526407	1
Caelo	Propylene glycol	1533405	2
Caelo	Propylene glycol	15335403	15
Caesar & Loretz GmbH Hi...	Propylene glycol	15335403	1
19.04.2016 Phoenix 2,84EUR	Propylene glycol	15335403	1
Beiersdorf AG	Propylene glycol	15335404	1
Caesar & Loretz GmbH Hi...	Propylene glycol	15335404	1
Caelo	Propylene glycol	15335404	15
Sanacorp	Propylene glycol	15335404	1
Sanacorp	Propylene glycol	15335405	1
lamotte/ Phoenix	Propylene glycol	15335405	2
Caelo	Propylene glycol	15335405	24
Noweda	Propylene glycol	15335405	1
Caesar&Loretz GmbH	Propylene glycol	15335405	1
phönix	Propylene glycol	15335405	1
Caelo	Propylene glycol	15335406	22
Caelo	Propylene glycol	15335407	16
Fagron	Propylene glycol	15335407	1
Caelo	Propylene glycol	15335408	12
Caelo	Propylene glycol	15355107	1
Caelo	Propylene glycol	15355108	10
Euro OTC	Propylene glycol	15355108	1
Caelo	Propylene glycol	153565101	1
Caelo	Propylene glycol	1540502	1
Caelo	Propylene glycol	15430502	12
Caelo	Propylene glycol	15413409	1
Caelo	Propylene glycol	1540513	1
phönix	Propylene glycol	15430502	1
Cealo/Phönix	Propylene glycol	15430502	1
sanacorp, 3,87EUR , 12.09...	Propylene glycol	154305025	1
Caelo	Propylene glycol	15430503	26
Caelo	Propylene glycol	15430504	12
Bombastus	Propylene glycol	15430504	2
phönix	Propylene glycol	15430504	2
Caelo	Propylene glycol	15430505	4
Caelo	Propylene glycol	15430506	8
Wolf-Apotheke	Propylene glycol	15430506	1
Caesar & Loretz GmbH Hi...	Propylene glycol	15430508	1
Caelo	Propylene glycol	15430508	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Propylene glycol	15430508	1
Caelo	Propylene glycol	15430512	3
Caelo	Propylene glycol	15430511	1
Lamotte / Noweda	Propylene glycol	15430512	2
Caelo	Propylene glycol	15430513	6
Alliance H.	Propylene glycol	15430513	1
Caelo	Propylene glycol	15430514	5
Caelo	Propylene glycol	15439502	1
Caelo	Propylene glycol	1562408	1
EuroPTC/Phönix	Propylene glycol	15901013-01	1
Fagron	Propylene glycol	15A07B04302970	1
Fagron	Propylene glycol	15A07B04	3
Fagron	Propylene glycol	15A07-B04	3
Fagron	Propylene glycol	15A07-B04-302969	7
Gehe	Propylene glycol	15A07-B04-302969	1
Fagron	Propylene glycol	15A07-B04-302970	18
Phönix Gotha	Propylene glycol	15A07-B04-302970	2
Caelo	Propylene glycol	15A07-B04-302970	3
Fagron	Propylene glycol	15B24-B40	3
Fagron	Propylene glycol	15B24-B40-307379	2
Caelo	Propylene glycol	15F25-B02-311890	1
Fagron	Propylene glycol	15F25-B02-311890	5
Fagron	Propylene glycol	15G28-B03	2
Fagron	Propylene glycol	15G28-B03-311889	3
Caelo	Propylene glycol	15G28-B03-311889	1
Fagron	Propylene glycol	15G28-B03-314768	19
Fagron	Propylene glycol	15L23-B01-3108036	1
Fagron	Propylene glycol	15G28-B03-314769	1
Fagron	Propylene glycol	15L23-B01-318036	12
Caelo	Propylene glycol	15L23-B01-318036	1
Euro OTC	Propylene glycol	16016014-01	1
Caelo	Propylene glycol	16021506	1
Euro OTC	Propylene glycol	1603034	1
Caelo	Propylene glycol	1603034-01	1
Euro OTC	Propylene glycol	1603034-01	16
Euro OTC	Propylene glycol	160303401	1
Euro OTC	Propylene glycol	160303401-01	1
Euro OTC	Propylene glycol	16060114-01	1
Euro OTC	Propylene glycol	1606014-01	9
Euro OTC	Propylene glycol	1610031	1
Euro OTC	Propylene glycol	1610031-01	10
Caelo	Propylene glycol	16105403	34
BASF/Gehe	Propylene glycol	1610031-01	1
Carlo/ Gehe	Propylene glycol	16105403	1
Calo/Sanacorp	Propylene glycol	16105403	1
Caesar & Loretz/Noweda	Propylene glycol	16105403	1
Caelo	Propylene glycol	16105404	5
Caelo	Propylene glycol	162251207	1
Caelo	Propylene glycol	1615403	1
Caelo	Propylene glycol	1625120	1
Caelo	Propylene glycol	16251201	27
caesar&loretz GmbH,Phön...	Propylene glycol	16251201	1
Caesar & Loretz GmbH	Propylene glycol	16251201	1
Caelo	Propylene glycol	16251202	25
phönix	Propylene glycol	16251202	4
Hedinger vdL	Propylene glycol	16251202	2
Gehe	Propylene glycol	16251202	2
Caelo	Propylene glycol	16251203	7
Noweda	Propylene glycol	16251202	1
Caelo	Propylene glycol	16251205	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix	Propylene glycol	16251205	1
Caelo	Propylene glycol	16251206	25
Caelo	Propylene glycol	16251207	6
Gehe	Propylene glycol	16251207	1
Fagron	Propylene glycol	16C10-B04	5
Caelo	Propylene glycol	16851205	1
Caelo	Propylene glycol	16551207	1
Caelo	Propylene glycol	16251208	2
Fagron	Propylene glycol	16C-10-B04	1
Fagron	Propylene glycol	16C10-B04-321081	19
Fagron	Propylene glycol	16c10-b04-321081	1
Fargon	Propylene glycol	16C10-B04-321081	1
Euro OTC	Propylene glycol	16C10-B04-321081	1
Caelo	Propylene glycol	16C10-B04-321659	1
Fagron	Propylene glycol	16C10-B04-323302	2
Caelo	Propylene glycol	16C10-B04-323302	2
Fagron	Propylene glycol	16C10-B04-321081	1
Fagron	Propylene glycol	16E09-B02	2
Fagron	Propylene glycol	16E09-B02-3024843	1
Fagron	Propylene glycol	16E09-B02-324	2
Fagron	Propylene glycol	16E09-B02-324843	4
Phönix	Propylene glycol	16E09-B02-324843	1
Fagron	Propylene glycol	16E09-B02-325331	1
Hubertus Apotheke	Propylene glycol	16G04-B01-328151	1
Fagron	Propylene glycol	16G04-B01-328151	3
Fagron	Propylene glycol	16G04B01328151	1
GEHE	Propylene glycol	16G04-B01-328151	1
Caelo	Propylene glycol	16G04B01328151	1
Fagron	Propylene glycol	16G04-B01-329687	1
Caelo	Propylene glycol	17041501	1
Caelo	Propylene glycol	17053801	17
Euro OTC	Propylene glycol	1703024-01	1
Caelo	Propylene glycol	170131PG/16105403	1
phönix	Propylene glycol	17053801	1
Caelo	Propylene glycol	17053802	3
Caelo	Propylene glycol	17053803	1
Caelo	Propylene glycol	17053804	14
GEHE	Propylene glycol	17053804	1
phönix	Propylene glycol	17053806	2
Caelo	Propylene glycol	17053806	19
Caesar und Loretz GmbH	Propylene glycol	17053806	2
Caeolo / Phönix	Propylene glycol	17053806	2
Caelo	Propylene glycol	17053807	3
Caelo	Propylene glycol	17053809	5
GEHE	Propylene glycol	17053809	2
Caelo	Propylene glycol	17053810	1
Caelo	Propylene glycol	17053814	2
Euro OTC	Propylene glycol	1708002-01	3
Caelo	Propylene glycol	170808PG/16251206	1
Euro OTC	Propylene glycol	1711023-01	2
Fagron	Propylene glycol	17082015A	1
Caelo	Propylene glycol	17275502	7
Caelo	Propylene glycol	171106PG/16251202	1
GEHE	Propylene glycol	17321703	1
Caelo	Propylene glycol	17321703	5
Caelo	Propylene glycol	17275504	1
Caelo	Propylene glycol	17321704	5
Caelo	Propylene glycol	1793451	3
Fagron	Propylene glycol	17A06-B03	1
Fagron	Propylene glycol	17A06-B03-332091	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Propylene glycol	17a06-b03-332091	1
GEHE	Propylene glycol	17A06-B03-332091	1
Fagron	Propylene glycol	17A06-B03-334896	1
Fagron	Propylene glycol	17A06-B03-332091	1
Fagron	Propylene glycol	17E22B10336830	1
Fagron	Propylene glycol	17H03-B04-340573	1
Fagron	Propylene glycol	17E22-B10-338768	1
GEHE	Propylene glycol	17H03-B04-340573	1
Caelo	Propylene glycol	180130PG/17275502	1
Fagron	Propylene glycol	18061321	2
Apotheke-Ansatz	Propylene glycol	19/15	1
Fagron	Propylene glycol	19111421	4
Fagron	Propylene glycol	18061311	1
Caelo	Propylene glycol	18064902	1
Caelo	Propylene glycol	18064901	1
Caelo	Propylene glycol	19111421	1
Euro OTC	Propylene glycol	200749	1
Fagron	Propylene glycol	200749-0001	1
Sonnen-Apotheke Wefensl...	Propylene glycol	20/16	1
Fagron	Propylene glycol	23011420	1
Caelo	Propylene glycol	26021320	2
Caelo	Propylene glycol	26041402	1
Fagron	Propylene glycol	270114A	1
Caelo	Propylene glycol	27091402	2
Lamotte/Phönix	Propylene glycol	2768	2
Caelo	Propylene glycol	28071503	1
Bombastus	Propylene glycol	289078	3
Euro OTC	Propylene glycol	289078	2
Fagron	Propylene glycol	280714C	1
Bombastus	Propylene glycol	286382	4
Sonnen-Apotheke Wefensl...	Propylene glycol		1
Caelo	Propylene glycol	29091403	1
Bombastus	Propylene glycol	291103	1
Bombastus	Propylene glycol	293947	8
Bombastus	Propylene glycol	295111	5
Bombastus	Propylene glycol	296424	4
Bombastus	Propylene glycol	297626	10
Bombastus	Propylene glycol	2983947	1
Bombastus	Propylene glycol	299106	14
Bombastus	Propylene glycol	301162	18
Caelo	Propylene glycol	300915	1
Bombastus	Propylene glycol	303788	19
Bombastuts/Gehe	Propylene glycol	307185	2
Caelo	Propylene glycol	307185	1
Phoenix	Propylene glycol	307185	1
Bombastus	Propylene glycol	307185	1
Bombastus	Propylene glycol	309030	3
Bombastus	Propylene glycol	309031	1
Bombastus	Propylene glycol	311416	1
Apotheke/Defektur Misch...	Propylene glycol	35/15	1
Sanacorp	Propylene glycol	50/AKD3	1
Sanacorp	Propylene glycol	49/AKD1	1
Caelo	Propylene glycol	31256708	1
Caelo	Propylene glycol	4248702	1
Caelo	Propylene glycol	5262402	1
Caelo	Propylene glycol	54/AKD7	1
Euro OTC	Propylene glycol	5315HH9041	2
Fagron	Propylene glycol	5115Q-02110	1
Henry Lamotte Oils GmbH	Propylene glycol	5471401	1
Henry Lamotte Oils GmbH...	Propylene glycol	5671401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils GmbH	Propylene glycol	5671401	10
Lamotte	Propylene glycol	5671401	10
Otto FischerGmbH / Alli...	Propylene glycol	5671401	1
lamotte	Propylene glycol	5671401	2
Henry Lamotte/Anzag	Propylene glycol	5671401	2
Caelo	Propylene glycol	5671401	7
Lamotte/Noweda	Propylene glycol	5671401	1
Lamotte / Alliance Heal...	Propylene glycol	5671401	1
Europ. AB	Propylene glycol	5671401	1
Sanacorp	Propylene glycol	5671401	1
Caelo	Propylene glycol	58/AKD16	1
Caelo	Propylene glycol	6509HH9041	1
Euro OTC	Propylene glycol	6509HH9041	1
Henry Lamotte/Alliance	Propylene glycol	8002520002	1
Henry Lamotte Oils	Propylene glycol	8002520002	1
AHD Lamotte GmbH	Propylene glycol	8002520002	1
Lamotte	Propylene glycol	8002520002	3
Henry Lamotte/Noweda	Propylene glycol	8002520002	1
Henry Lamotte / Noweda	Propylene glycol	8002520002	1
Henry Lamotte Oils GmbH...	Propylene glycol	8002520002	1
Lamotte/Sanacorp	Propylene glycol	8002520002	1
Alliance/Henry	Propylene glycol	8002520002	1
gehe	Propylene glycol	8002520002	1
Henry Lamotte	Propylene glycol	8002520002	1
Caelo	Propylene glycol	8002520002	1
henry lamotte/gehe	Propylene glycol	8002520002	1
Henry Lamotte	Propylene glycol	8003093001	3
Lamotte	Propylene glycol	8003093001	2
Lamotte GmbH/ Kehr	Propylene glycol	8003093001	1
Sanacorp	Propylene glycol	AKD70	1
Caelo	Propylene glycol	83/2014	1
Caelo	Propylene glycol		2
Sanacorp	Propylene glycol		1
Bombastus	Propylene glycol	268295	1
Euro OTC	Propylene glycol	q603034-01	1
ACM	Propylene glycol	KNC0776J15	1
Fagron	Propylene glycol	TQ131410	2
Bombastus	Propylene glycol	265358	1
Bombastus	Propylene glycol	270359	1
Bombastus	Propylene glycol	7/14	1
Caelo	Propylene glycol	113Q-02110	1
Lamotte 004	Propylene glycol	122882	1
Noweda / Henry Lamotte	Propylene glycol	122882	1
Celo	Propylene glycol	11607-2	1
Henry Lamotte,Phönix	Propylene glycol	123095	2
Großhandel	Propylene glycol	122882	1
Lamotte	Propylene glycol	122690	1
Lamotte	Propylene glycol	123095	1
Henry Lamotte	Propylene glycol	123095	4
Henry Lamotte Oils GmbH	Propylene glycol	123095	2
Bombastus	Propylene glycol	123095	1
Henry Lamotte oils GmbH	Propylene glycol	123095	1
Henry Lam./ Phönix	Propylene glycol	123095	1
henry lamotte	Propylene glycol	125099	1
Caelo	Propylene glycol	2101302	1
Hedinger	Propylene glycol	29595B	1
Hedinger	Propylene glycol	29625B	1
Caelo	Propylene glycol	5061201	1
Caelo	Propylene glycol	6021301	1
Caelo	Propylene glycol	611211	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils	Propylene glycol	7021401	1
Caelo	Propylene glycol	70214PJK02	1
Fagron	Propylene glycol	8121120	1
Fagron	Propylene glycol	9031321	1
Fagron	Propylene glycol	9051220	7
Phönix Gotha	Propylene glycol	10210312	1
henry lamotte	Propylene glycol	1000001132	1
Caelo	Propylene glycol	10210307	1
Alliance Healthcare	Propylene glycol	10274303	1
Caelo	Propylene glycol	10274305	2
Caelo	Propylene glycol	10371307	1
Caelo	Propylene glycol	11021404	1
Euro OTC	Propylene glycol	1109015-01	1
Noweda	Propylene glycol	1112004-01	2
Caelo	Propylene glycol	11125804	1
Caelo	Propylene glycol	11300502	2
Caelo	Propylene glycol	11125811	1
Caelo	Propylene glycol	11125812	1
Caelo	Propylene glycol	11300505	2
Gehe	Propylene glycol	11300506	1
Caelo	Propylene glycol	11300506	5
Caesar & Loretz GmbH, 4...	Propylene glycol	11300505	1
Kehr	Propylene glycol	11300506	1
Noweda	Propylene glycol	11300506	1
Phönix	Propylene glycol	1130505	1
Euro OTC	Propylene glycol	1136HH9041	1
Caelo	Propylene glycol	11372005	3
Noweda	Propylene glycol	11372005	2
Erika apotheke	Propylene glycol	11372006	1
Anzag ek: 3.06EUR 15.08.12	Propylene glycol	11372006	3
noweda	Propylene glycol	11372006	1
Caelo	Propylene glycol	11372006	12
Caesar & Loretz GmbH	Propylene glycol	11372006	1
Caelo	Propylene glycol	11372007	8
Noweda	Propylene glycol	11372007	1
Euro OTC	Propylene glycol	11372007	2
Gehe	Propylene glycol	11372007	3
Caesar & Loretz GmbH, 4...	Propylene glycol	11372007	1
Phönix 28.08.2012	Propylene glycol	11372007	1
Caelo	Propylene glycol	11372009	5
Noweda	Propylene glycol	11372009	1
Gehe	Propylene glycol	11372009	1
Caelo	Propylene glycol	11372010	6
Phönix Gotha	Propylene glycol	11372011	1
Caelo	Propylene glycol	11372011	4
Euro OTC	Propylene glycol	11C07-N06	1
Fagron	Propylene glycol	11F10-N04	1
v.d.l	Propylene glycol	11h04-n01	1
Fagron	Propylene glycol	11H04-N01	3
Caelo	Propylene glycol	11H04-N01	1
Fagron	Propylene glycol	11K28-N01	10
Caelo	Propylene glycol	11K28-N01	1
Fagron	Propylene glycol	11L22-N14	7
Caelo	Propylene glycol	1206038-01	1
Euro OTC	Propylene glycol	1206038-01	4
Sanacorp	Propylene glycol	1206038-01	1
Euro OTC	Propylene glycol	120603801	1
Noweda	Propylene glycol	12088205	3
Caelo	Propylene glycol	12088205	3
Caelo	Propylene glycol	12088206	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Propylene glycol	12088205	1
Herbasinica	Propylene glycol	12088206	2
Caelo	Propylene glycol	12088207	16
Noweda	Propylene glycol	12088207	2
Phytocomm	Propylene glycol	12088207	1
Eigenherstellung	Propylene glycol	120912PML01	1
Caelo	Propylene glycol	1208827	1
Fagron	Propylene glycol	120C004-277798	1
Caelo	Propylene glycol	12100907	1
Caelo	Propylene glycol	12111304	1
Caesar & Loretz GmbH, 4...	Propylene glycol	12169802	5
Caelo	Propylene glycol	12169802	17
Anzag 25.01.2013	Propylene glycol	12169802	2
Noweda	Propylene glycol	12169802	2
Gehe	Propylene glycol	12169802	1
Caelo	Propylene glycol	12169804	8
Caelo	Propylene glycol	12372403	7
Alliance Healthcare	Propylene glycol	12372403	1
Caelo	Propylene glycol	12372402	2
Fagron	Propylene glycol	12372403	1
Caelo	Propylene glycol	12372404	14
Anzag 19.03.13 ek: 3,03EUR	Propylene glycol	12372404	2
Caesar & Loretz GmbH	Propylene glycol	12372404	1
Caelo	Propylene glycol	12372405	16
Gehe	Propylene glycol	12372405	2
Phönix	Propylene glycol	12372405	1
Gehe 00.07.2013	Propylene glycol	12372405	1
Caelo	Propylene glycol	1288205	1
Caelo	Propylene glycol	1272404	1
Fagron	Propylene glycol	12C06-N04	6
Caelo	Propylene glycol	12C06-N04	1
Fagron	Propylene glycol	12G13-N07	31
Fagron	Propylene glycol	12g13-n07	3
Caelo	Propylene glycol	12G13-N07	4
Fagron	Propylene glycol	12G13-N07	5
Fagron	Propylene glycol	12G13-No7	1
Gehe	Propylene glycol	12J04-N06	2
Fagron	Propylene glycol	12K13-N05	15
Gehe	Propylene glycol	12K13-N05	2
Fagron	Propylene glycol	12k13-n05	2
Caelo	Propylene glycol	12K13-N05	1
Fargon	Propylene glycol	12K13-N05	2
Fagron	Propylene glycol	1301C003	3
Fagron	Propylene glycol	1301C003-277799	1
Fagron	Propylene glycol	1301C004	5
Fagron	Propylene glycol	1301C004-277798	16
Caelo	Propylene glycol	1301C004-277798	2
Phönix Gotha	Propylene glycol	1301C004-277798	1
EuRho	Propylene glycol	1302002-01	1
Phoenix	Propylene glycol	1301C004-277798	1
Euro OTC	Propylene glycol	1302002-01	11
Noweda	Propylene glycol	1302002-01	1
Euro OTC	Propylene glycol	130200201	2
Fagron	Propylene glycol	1304C009-28	1
Caelo	Propylene glycol	1304C009-281100	2
Fagron	Propylene glycol	13040010-281099	1
Fagron	Propylene glycol	1304C010	11
Fagron	Propylene glycol	1304C-010	1
Gehe	Propylene glycol	1304C010-281099	2
Fagron	Propylene glycol	1304C010-281099	21

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Propylene glycol	1304c010-281099	4
Fagron	Propylene glycol	1304C010281099	1
Phoenix	Propylene glycol	1304C010-281099	1
Fagron	Propylene glycol	1304C010	1
Caelo	Propylene glycol	13071303	1
Caelo	Propylene glycol	13100906	12
Euro OTC	Propylene glycol	130702101	1
Caesar/Noweda	Propylene glycol	13100906	1
Caelo	Propylene glycol	13100907	23
AHD	Propylene glycol	13100907	1
Phoenix	Propylene glycol	13100907	1
Cealo	Propylene glycol	13100907	2
Caelo	Propylene glycol	13100908	37
Gehe	Propylene glycol	13100908	3
Caelo	Propylene glycol	1316003	1
Caelo	Propylene glycol	13167003	26
Noweda	Propylene glycol	1313HH9041	1
Caelo	Propylene glycol	1310098	1
Caesar /phönix	Propylene glycol	13167003	1
Noweda	Propylene glycol	13167003	1
canea, anzag	Propylene glycol	13167003	3
Phönix	Propylene glycol	13167003	2
Gehe	Propylene glycol	13167003	2
Caelo	Propylene glycol	13167004	28
Phönix	Propylene glycol	13167004	1
Caelo	Propylene glycol	13167005	1
Caelo	Propylene glycol	13167006	7
Caesar/Noweda	Propylene glycol	13167006	2
Phönix	Propylene glycol	13167006	1
Sanacorp	Propylene glycol	13167006	1
Caelo	Propylene glycol	1316Q-02110	1
Noweda	Propylene glycol	13256704	1
Gehe	Propylene glycol	13256704	2
ANZAG 20.01.2014	Propylene glycol	13256704	1
Phoenix	Propylene glycol	13256705	1
Fagron	Propylene glycol	13E22B03	1
Fagron	Propylene glycol	13E22B0328345	1
gehe	Propylene glycol	13E22-B03-283450	1
Anzag	Propylene glycol	13E22-B03-283450	1
Fagron	Propylene glycol	13E22-B03-283451	2
Euro OTC	Propylene glycol	140101001	1
Fagron	Propylene glycol	14051802	2
Phönix	Propylene glycol	14248706	1
Phönix/	Propylene glycol	14248707	3
Caesar&lloretz,Phoenix	Propylene glycol	15069202	1
Caelo	Propylene glycol	150771	1
Phoenix	Propylene glycol	18061321	1
Fagron	Propylene glycol	18081121	2
Caelo	Propylene glycol	15A07-B04-302969	1
Caelo	Propylene glycol	18041302	1
Caelo	Propylene glycol	1804212110	1
Fischbar	Propylene glycol	18091201	1
Gehe	Propylene glycol	18081121	1
Fagron	Propylene glycol	181013prop	1
Caelo	Propylene glycol	19041307	1
Caelo	Propylene glycol	20071201	1
Fagron	Propylene glycol	200749-B04-289597	1
Caelo	Propylene glycol	20111203	1
Fagron	Propylene glycol	20122013B	1
Caelo	Propylene glycol	20140612-6	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Propylene glycol	24041221	2
Fagron	Propylene glycol	2112212110	1
Fagron	Propylene glycol	24041221	9
Gehe	Propylene glycol	24041221	1
Caelo	Propylene glycol	24121301	1
Caelo	Propylene glycol	25011302	1
Caelo	Propylene glycol	25031301	1
Fagron	Propylene glycol	26021320	6
Caelo	Propylene glycol	26061302	1
Alliance Healthcare	Propylene glycol	272684	1
Bombastus	Propylene glycol	277548	4
Caelo	Propylene glycol	28021301	1
Bombastus	Propylene glycol	281569	3
Bombastus	Propylene glycol	283937	1
Caelo	Propylene glycol	31051317	1
Caelo	Propylene glycol	42488055	1
Caelo	Propylene glycol	2913M-02110	1
Fagron	Propylene glycol	3959/V02/15	1
Caelo	Propylene glycol	3512M-02110	1
Fagron	Propylene glycol		2
phönix	Propylene glycol	AKD30	1
Henry Lamotte	Propylene glycol	5671401	1
Fagron	Propylene glycol	Ch12J04-N06	1
Caelo	Propylene glycol	Ch13100908	1
Phönix	Propylene glycol	Retest03083	1
Gehe	Propylene glycol	Test2	1
Henry Lamotte Oils	Propylene glycol	122880	1
Caelo	Propylene glycol	122882	1
Phagron	Propylene glycol		1
Henry Lamotte Oils	Propylene glycol	123095	3
Phönix	Propylene glycol	29595B	1
A.H.D.	Propylene glycol	1112004-01	1
Caelo	Propylene glycol	11300504	1
GEHE	Propylene glycol	11372007	1
Fagron	Propylene glycol	11h04-n01	1
Anzag	Propylene glycol	12G13N07	1
Gehe	Propylene glycol	12088207	1
Caelo	Propylene glycol	1301C003-277799	1
Fagron	Propylene glycol	1301c004	2
Fagron	Propylene glycol	1301C004-0005	1
Sanacorp	Propylene glycol	1302002-01	1
Fagron	Propylene glycol	1304c009-281100	1
Caelo	Propylene glycol	13100905	2
Fagron	Propylene glycol	13G09-B01-283452	1
Fagron	Propylene glycol	o903132121	1

- 88 609 spectra from 501 *Apo-Ident* customers from a total of 26 562 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Propylene glycol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all

relevant spectra of the various substances were compared with *Propylene glycol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	285	0	35 259
Type B	0	315	0	93 559
Type C	0	2716	27	88 609

The substance/substance group *Propylene glycol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.8947 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.0952 %)
Type C	100.0000 % (> 99.6586 %)	99.0157 % (> 98.9063 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30872	31663	3.34	47.77
30920	30920	0.00	49.63
31663	31663	0.00	50.90
31738	31738	0.00	48.33
32074	31663	5.31	54.29
32486	31663	6.75	57.29
32728	31663	5.63	48.51

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Psorisol base cream Audor[®]**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31109-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Psorisol base cream Audor[®]

Special notes

When selecting the *Psorisol base cream Audor[®]* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Psorisol base cream Audor[®]*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Audor	Psorisol base cr...	AP101943PB	31109	40	1404015
Audor	Psorisol base cr...	AP45001.2PB	31674	40	20141230
Audor	Psorisol base cr...	AP55201PB	32275	40	20160204

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Psorisol base cream Audor®*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Psorisol base cream Audor®*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Audor Pharma	Psorisol base cream Audor®	AP102640PB	30281	40
Audor Pharma	Psorisol base cream Audor®	AP101943PB	31233	40
Audor Pharma	Psorisol base cream Audor®	AP45001.2PB	31674 [†]	20
Audor Pharma	Psorisol base cream Audor®	AP45001.2PB	32274	40
Audor Pharma	Psorisol base cream Audor®	AP70601PB	34037	40

- 93 694 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 25 spectra from 10 *Apo-Ident* customers from 13 batches from the substance/substance group *Psorisol base cream Audor®*.
- These include spectra of independent samples from 9 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Audor Pharma	Psorisol base cream Audor®	101943	1
Audor Pharma	Psorisol base cream Audor®		2
Audor Pharma	Psorisol base cream Audor®	AP102844PB	3
Audor Pharma	Psorisol base cream Audor®	102844	2
Audor Pharma	Psorisol base cream Audor®	AP55201PB	6
Audor Pharma	Psorisol base cream Audor®	AP101943PB	3
Audor Pharma	Psorisol base cream Audor®	102891PB	1
Audor Pharma	Psorisol base cream Audor®	102891	1
Audor Pharma	Psorisol base cream Audor®	AP100472PB	2
Audor Pharma	Psorisol base cream Audor®	AP11943PB	1
Audor Pharma	Psorisol base cream Audor®	AP102891PB	1
Anzag	Psorisol base cream Audor®	AP102891PB	2

- 91 327 spectra from 501 *Apo-Ident* customers from a total of 26 996 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Psorisol base cream Audor®* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Psorisol base cream Audor®* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	180	0	93 694
Type C	0	18	7	91 327

The substance/substance group *Psorisol base cream Audor®* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.6667 %)
Type C	100.0000 % (> 99.6589 %)	72.0000 % (> 60.0000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several

new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31109	31109	0.00	22.63
31674	31109	4.67	21.47
32275	31109	6.50	22.19

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Rose oil (genuine / artificial)**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30111-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Rose oil (genuine / artificial); Oleum rosae (verum / artificiale)

Special notes

When selecting the *Rose oil (genuine / artificial)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Rose oil (genuine / artificial)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Rose oil (genuin...	12164308	31021	40	1502254
Caelo	Rose oil (genuin...	11287117	31022	40	20111011
Fagron	Rose oil (genuin...	13J17-B01	31290	40	1507614
Caelo	Rose oil (genuin...	13335012	31301	20	20131009
Caelo	Rose oil (genuin...	13252710	31412	40	1502252

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Rose oil (genuin...	13335012	31425	40	20131009
Fagron	Rose oil (genuin...	14I25-B01-302286	31752	40	20141006
Fagron	Rose oil (genuin...	14I25-B01-302286	31842	40	20141006

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 300 spectra of 8 reference samples from the substance/substance group *Rose oil (genuine / artificial)*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 6 different batches.
- 35 244 spectra from a total of 834 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 243 spectra from 7 reference samples of the substance/substance group *Rose oil (genuine / artificial)*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Rose oil (genuine / artifici...	09185909	30209	43
Caelo	Rose oil (genuine / artifici...	13335012	31425 [†]	20
Caelo	Rose oil (genuine / artifici...	13335012	31426	60
Fagron	Rose oil (genuine / artifici...	14I25-B01-302286	31752 [†]	20
Fagron	Rose oil (genuine / artifici...	14I25-B01-302286	31842 [†]	20
Caelo	Rose oil (genuine / artifici...	162311	32332	40
Caelo	Rose oil (genuine / artifici...	173537	33907	40

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*).

- 93 631 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 127 spectra from 67 *Apo-Ident* customers from 58 batches from the substance/substance group *Rose oil (genuine / artificial)*.
- These include spectra of independent samples from 55 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Rose oil (genuine / artifici...	11287117	2
Caelo	Rose oil (genuine / artifici...	11287114	1
Fagron	Rose oil (genuine / artifici...	12C19-N05	3
Fagron	Rose oil (genuine / artifici...	7091301	1
Fagron	Rose oil (genuine / artifici...	12H16-B01-286704	3
Fagron	Rose oil (genuine / artifici...	12G25-N02	6
Fagron	Rose oil (genuine / artifici...	13C21-N04	10
Phönix	Rose oil (genuine / artifici...	13C21-N04	2
Fagron	Rose oil (genuine / artifici...	698884	3
Caelo	Rose oil (genuine / artifici...	13252710	3
Fagron	Rose oil (genuine / artifici...	12H16-B01-286704	3
Caelo	Rose oil (genuine / artifici...	13252723	1
Caelo	Rose oil (genuine / artifici...	14300712	4
Fagron	Rose oil (genuine / artifici...	14I25-B01	3
Caelo	Rose oil (genuine / artifici...	15197202	11
Fagron	Rose oil (genuine / artifici...	12G25-N01	1
Caelo	Rose oil (genuine / artifici...	11287123	1
Fagron	Rose oil (genuine / artifici...	12H16-B01	3
Caelo	Rose oil (genuine / artifici...	13252701	2
Caelo	Rose oil (genuine / artifici...	13252722	4
Caelo	Rose oil (genuine / artifici...	1352710	1
Caelo	Rose oil (genuine / artifici...	13335006	1
Fagron	Rose oil (genuine / artifici...	13J17-801	1
Fagron	Rose oil (genuine / artifici...	13J17-B01	3
Fagron	Rose oil (genuine / artifici...	13C21-N074	1
Fagron	Rose oil (genuine / artifici...	13J17-B01-293574	3
Fagron	Rose oil (genuine / artifici...	14/25-B01-302286	1
Fagron	Rose oil (genuine / artifici...	14/25-B01-309675	1
Fagron	Rose oil (genuine / artifici...	13J17-B01-296955	1
Fagron	Rose oil (genuine / artifici...	14I25-B01	1
Caelo	Rose oil (genuine / artifici...	14300703	1
Fagron	Rose oil (genuine / artifici...	1425-B01-302287	1
Fagron	Rose oil (genuine / artifici...	14I25-B01-309676	1
Caelo	Rose oil (genuine / artifici...	14300719	1
Fagron	Rose oil (genuine / artifici...	15/25-B01-302286	1
Caelo	Rose oil (genuine / artifici...	15187202	1

continued on the next page

Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Rose oil (genuine / artifici...	17101412	1
Caelo	Rose oil (genuine / artifici...	19011502	1
Fagron	Rose oil (genuine / artifici...	2188	1
Fagron	Rose oil (genuine / artifici...	713371	1
Fagron	Rose oil (genuine / artifici...	1425-01-309676	1
Caelo	Rose oil (genuine / artifici...	15197206	1
Caelo	Rose oil (genuine / artifici...	15197211	8
Fagron	Rose oil (genuine / artifici...	15G29-B12	1
Frey & Lau GmbH/Noweda	Rose oil (genuine / artifici...	15G29-B12	1
Caelo	Rose oil (genuine / artifici...	15G29-B12-312961	1
Caelo	Rose oil (genuine / artifici...	15G29-B12-314833	1
Fagron	Rose oil (genuine / artifici...	15G29-B12-321984	3
Fagron	Rose oil (genuine / artifici...	15G29-B12-327263	1
Fagron	Rose oil (genuine / artifici...	15G29-B12-314833	1
Caelo	Rose oil (genuine / artifici...	16231111	1
Fagron	Rose oil (genuine / artifici...	15G29-B12-327364	1
Caelo	Rose oil (genuine / artifici...	17353703	1
Caelo	Rose oil (genuine / artifici...	16231113	1
Caelo	Rose oil (genuine / artifici...	11178601	1
Caelo	Rose oil (genuine / artifici...	13044608	2
Fagron	Rose oil (genuine / artifici...	11E04-N04	1
Caelo	Rose oil (genuine / artifici...	12164306	1
Fagron	Rose oil (genuine / artifici...	11H05-N02	1
FARGON	Rose oil (genuine / artifici...	12c19-n05	1
Fagron	Rose oil (genuine / artifici...	511I-02011	1
Phönix	Rose oil (genuine / artifici...	12C19-N05	1
Phönix	Rose oil (genuine / artifici...	12G25-N02	2
Gehe 29.06.2013	Rose oil (genuine / artifici...	12G25-N02	1
Fagron	Rose oil (genuine / artifici...	887963	1

- 91 225 spectra from 501 *Apo-Ident* customers from a total of 26 951 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Rose oil (genuine / artificial)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Rose oil (genuine / artificial)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	300	0	35 244
Type B	0	243	0	93 631
Type C	0	121	6	91 225

The substance/substance group *Rose oil (genuine / artificial)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.5309 %)
Type C	100.0000 % (> 99.6586 %)	95.2756 % (> 92.9134 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31021	31021	0.00	84.87
31022	31412	9.63	63.69
31290	31290	0.00	143.65
31301	31021	9.72	78.43
31412	31412	0.00	59.68
31425	31021	5.71	80.34
31752	31290	3.52	140.67
31842	31290	3.70	138.37

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Sage oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30578-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Sage oil; Oleum salviae; Salvia officinalis oleum

Special notes

When selecting the *Sage oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Sage oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sage oil	13146901	31597	40	20130522
Bombastus	Sage oil	305303	33292	40	1709470
Caelo	Sage oil	16051006	33601	40	1803004

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Sage oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 180 spectra from 5 reference samples of the substance/substance group *Sage oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Sage oil	10248202	30170	40
Caelo	Sage oil	11287302	30578	40
Caelo	Sage oil	13146901	31597 [†]	20
Bombastus	Sage oil	305303	33453	40
Caelo	Sage oil	16051012	33949	40

- 93 694 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 99 spectra from 57 *Apo-Ident* customers from 42 batches from the substance/substance group *Sage oil*.
- These include spectra of independent samples from 41 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Sage oil	13146903	2
Caelo	Sage oil	13146901	8
Caelo	Sage oil	13146907	1
Caelo	Sage oil	13146914	6
Caelo	Sage oil	14345512	1
Caelo	Sage oil	13146908	2
Caelo	Sage oil	14345517	8
Caelo	Sage oil	17310001	1
Caelo	Sage oil	16051012	1
Caelo	Sage oil	10248206	2
Caelo	Sage oil	10248215	1
Anzag 10.8.12 ek:4,42	Sage oil	11287305	1
Anzag, 11.09.2012	Sage oil	11287305	1
Caelo	Sage oil	11287306	1
Caelo	Sage oil	11287307	1
Caelo	Sage oil	11287312	1
Caelo	Sage oil	11287314	2
Caelo	Sage oil	11287311	2
Caelo	Sage oil	11287315	1
Caesar & Loretz GmbH, 4...	Sage oil	11287317	1
Caelo	Sage oil	11287317	11
Caelo	Sage oil	11287318	1
Gehe	Sage oil	11287318	1
Caelo	Sage oil	11287320	3
Caelo	Sage oil	11287323	7
Caelo	Sage oil	13146904	1
gehe 19.11.2013	Sage oil	13146901	2
Caelo	Sage oil	21111301	1
Caelo	Sage oil	28051312	1
Sanacorp	Sage oil	282545	1
Caelo	Sage oil	31011401	1
Taoasis	Sage oil	498780-106211	1
Toasis	Sage oil	526975-110409	1
Taoasis GmbH	Sage oil	526975-110409	1
Taoasis	Sage oil	526975-111295	1
Caelo	Sage oil	5x11287317	1
Taoasis	Sage oil	530771-108279	1
Caelo	Sage oil	13146918	1
Bombastus	Sage oil	298189	1
Bombastus	Sage oil	274447	1
Caelo	Sage oil	1024214	1
Caelo	Sage oil	14345505	7
Caelo	Sage oil	14345511	3
Cealo	Sage oil	14345517	1
Caelo	Sage oil	1435517	1
Caelo	Sage oil	1702080S/14345517	1
Caelo	Sage oil	1704030S/14345517	1
Caelo	Sage oil	16051009	1

- 91 253 spectra from 501 *Apo-Ident* customers from a total of 26 967 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Sage oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Sage oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	169	11	93 694
Type C	0	86	13	91 253

The substance/substance group *Sage oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	93.8889 % (> 92.2222 %)
Type C	100.0000 % (> 99.6586 %)	86.8687 % (> 83.8384 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31597	33601	5.14	22.30
33292	33292	0.00	20.82
33601	33601	0.00	24.40

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Salicylic acid 50% (paraffin)
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30496-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Salicylic acid 50% (paraffin); Acidum salicylicum cum vaselino alba 50%

Special notes

When selecting the *Salicylic acid 50% (paraffin)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Salicylic acid 50% (paraffin)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Salicylic acid 5...	12305902	30883	40	20121107
Henry Lamotte	Salicylic acid 5...	5375300	30927	40	20121003
Caelo	Salicylic acid 5...	15231801	31727	40	AR-15-FG-015411-01
Caelo	Salicylic acid 5...	153729	31901	40	20151029*
Caelo	Salicylic acid 5...	153730	31902	40	20151029*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Salicylic acid 5...	153731	31903	40	20151103*
Henry Lamotte	Salicylic acid 5...	6001023001	32176	40	20150630
Caelo	Salicylic acid 5...	162761	32520	40	20160923*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 320 spectra of 8 reference samples from the substance/substance group *Salicylic acid 50% (paraffin)*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 8 different batches.
- 35 224 spectra from a total of 832 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 480 spectra from 13 reference samples of the substance/substance group *Salicylic acid 50% (paraffin)*.
- These include spectra of independent samples from 9 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
apomix	Salicylic acid 50% (paraffin)	04B270911	30363	20
apomix	Salicylic acid 50% (paraffin)	01A130212	30496	40
Caelo	Salicylic acid 50% (paraffin)	14107901	31357	60
Caelo	Salicylic acid 50% (paraffin)	15231801	31727 [†]	20
Caelo	Salicylic acid 50% (paraffin)	153729	31901 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Salicylic acid 50% (paraffin)	153730	31902 [†]	20
Caelo	Salicylic acid 50% (paraffin)	153731	31903 [†]	20
Caelo	Salicylic acid 50% (paraffin)	153732	31904	60
Caelo	Salicylic acid 50% (paraffin)	153733	31905	60
Caelo	Salicylic acid 50% (paraffin)	171454	33357	40
Caelo	Salicylic acid 50% (paraffin)	171055	33358	40
Caelo	Salicylic acid 50% (paraffin)	171727	33359	40
Caelo	Salicylic acid 50% (paraffin)	182109	34266	40

- 93 394 spectra from a total of 2154 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 686 spectra from 215 *Apo-Ident* customers from 164 batches from the substance/substance group *Salicylic acid 50% (paraffin)*.
- These include spectra of independent samples from 159 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Salicylic acid 50% (paraffin)	15231601	11
Caelo	Salicylic acid 50% (paraffin)	15231801	15
Caelo	Salicylic acid 50% (paraffin)	15247101	24
Caelo	Salicylic acid 50% (paraffin)	16121401	10
Caelo	Salicylic acid 50% (paraffin)	15305801	9
Caelo	Salicylic acid 50% (paraffin)	15305901	12
Caelo	Salicylic acid 50% (paraffin)	15373301	11
Caelo	Salicylic acid 50% (paraffin)	16276101	17
Henry Lamotte / Noweda	Salicylic acid 50% (paraffin)	6001023002	1
Caelo	Salicylic acid 50% (paraffin)	12201701	1
HENRY LAMOTTE OILS GmbH	Salicylic acid 50% (paraffin)	5375300	2
Phönix	Salicylic acid 50% (paraffin)	12201801	1
A.H.D.	Salicylic acid 50% (paraffin)	1228301	1
Caelo	Salicylic acid 50% (paraffin)	13160002	2
Caelo	Salicylic acid 50% (paraffin)	13087201	2
apomix	Salicylic acid 50% (paraffin)	1A150116	2
Apomix / 11	Salicylic acid 50% (paraffin)	1A180214	2
GEHE	Salicylic acid 50% (paraffin)	2A270217	1
APOMIX	Salicylic acid 50% (paraffin)	2A280116	2
GEHE	Salicylic acid 50% (paraffin)	3A070417	1
Noweda/apomix	Salicylic acid 50% (paraffin)	5A020615	1
Apomix / 11	Salicylic acid 50% (paraffin)	4A080515	1
Apomix	Salicylic acid 50% (paraffin)	6A150615	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix / Phönix	Salicylic acid 50% (paraffin)	6B110714	1
apomix	Salicylic acid 50% (paraffin)	5A070416	1
apomix	Salicylic acid 50% (paraffin)	6A270416	1
apomix	Salicylic acid 50% (paraffin)	7A230615	1
Noweda/Apomix	Salicylic acid 50% (paraffin)	8A250914	1
Caelo	Salicylic acid 50% (paraffin)		1
apomix	Salicylic acid 50% (paraffin)	9A140716	1
Noweda/Apomix	Salicylic acid 50% (paraffin)	7A230615	1
apomix	Salicylic acid 50% (paraffin)	9A150715	1
PKH Halle/Noweda	Salicylic acid 50% (paraffin)	9A171213	1
Apomix	Salicylic acid 50% (paraffin)	10A290716	1
Apomix/Gehe	Salicylic acid 50% (paraffin)	11A290716	1
Apomix / 11	Salicylic acid 50% (paraffin)	10A181114	1
apomix	Salicylic acid 50% (paraffin)	11A290716	1
apomix	Salicylic acid 50% (paraffin)	12A050916	1
Apomix/AHD	Salicylic acid 50% (paraffin)	12A050916	1
apomix	Salicylic acid 50% (paraffin)	12A130815	1
Caelo	Salicylic acid 50% (paraffin)	13278501	11
Caelo	Salicylic acid 50% (paraffin)	13278801	3
Caelo	Salicylic acid 50% (paraffin)	13278701	12
apomix	Salicylic acid 50% (paraffin)	13A121016	1
apomix	Salicylic acid 50% (paraffin)	13A240915	1
Caelo	Salicylic acid 50% (paraffin)	14017501	15
Caelo	Salicylic acid 50% (paraffin)	14017501	3
Caelo	Salicylic acid 50% (paraffin)	14107901	20
AHD 04.07.14 EK 5,91	Salicylic acid 50% (paraffin)	14107901	1
Caelo	Salicylic acid 50% (paraffin)	14108001	21
Raben-Apotheke Mildena	Salicylic acid 50% (paraffin)	14108201	1
Caelo	Salicylic acid 50% (paraffin)	14108201	3
Caelo	Salicylic acid 50% (paraffin)	14108101	3
Caelo	Salicylic acid 50% (paraffin)	14108301	4
Caelo	Salicylic acid 50% (paraffin)	14108402	1
Caelo	Salicylic acid 50% (paraffin)	14205101	10
Caelo	Salicylic acid 50% (paraffin)	14204302	1
Caelo	Salicylic acid 50% (paraffin)	14205001	3
Caelo	Salicylic acid 50% (paraffin)	14337201	10
Caelo	Salicylic acid 50% (paraffin)	14337301	13
Caelo	Salicylic acid 50% (paraffin)	14354101	3
Apomix	Salicylic acid 50% (paraffin)	14A031115	1
apomix	Salicylic acid 50% (paraffin)	14A121016	1
Apomix PKH	Salicylic acid 50% (paraffin)	15A251116	1
Caelo	Salicylic acid 50% (paraffin)	15372901	9
Caelo	Salicylic acid 50% (paraffin)	16197202	1
Caelo	Salicylic acid 50% (paraffin)	16197401	16
Hochhaus-Apotheke	Salicylic acid 50% (paraffin)	16197302	1
Caelo	Salicylic acid 50% (paraffin)	17017301	8
Lamotte/Gehe	Salicylic acid 50% (paraffin)	5763800	1
Brückenapotheke	Salicylic acid 50% (paraffin)	1715Q-07018	1
Henry Lamotte	Salicylic acid 50% (paraffin)	5763800	3
Caelo	Salicylic acid 50% (paraffin)	17145502	2
Caelo	Salicylic acid 50% (paraffin)	17326201	1
Lamotte	Salicylic acid 50% (paraffin)	6001862002	5
Henry Lamotte Oils GmbH...	Salicylic acid 50% (paraffin)	6004854001	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	6004854001	1
Phoenix (PKH Halle)	Salicylic acid 50% (paraffin)	2A040313	1
Apomix	Salicylic acid 50% (paraffin)	5A250713	1
Caelo	Salicylic acid 50% (paraffin)	13160101	5
apomix, PKH GmbH	Salicylic acid 50% (paraffin)	7A160913	1
Hersteller: apomix	Salicylic acid 50% (paraffin)	6A301111	1
Caelo	Salicylic acid 50% (paraffin)	13160102	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Salicylic acid 50% (paraffin)	13067201	1
Caelo	Salicylic acid 50% (paraffin)	12201901	5
Apomix/N/PZN:6639275	Salicylic acid 50% (paraffin)	9A171213	1
Caelo	Salicylic acid 50% (paraffin)	13087301	10
cealo	Salicylic acid 50% (paraffin)	13278501	3
Caelo	Salicylic acid 50% (paraffin)	13278601	10
Caelo	Salicylic acid 50% (paraffin)	13087401	7
Caelo	Salicylic acid 50% (paraffin)	14204902	1
Lamotte	Salicylic acid 50% (paraffin)	5065800	1
Henry Lamotte	Salicylic acid 50% (paraffin)	5375300	5
Caelo	Salicylic acid 50% (paraffin)	5375300	1
Henry Lamotte	Salicylic acid 50% (paraffin)	5448100	1
Lamotte/VDL	Salicylic acid 50% (paraffin)	5375300	1
Caelo	Salicylic acid 50% (paraffin)	s-100513-1	1
apomix	Salicylic acid 50% (paraffin)	2A260115	2
PKH / 11	Salicylic acid 50% (paraffin)	2A280116	1
Gehe	Salicylic acid 50% (paraffin)	14017501	1
Phönix Gotha	Salicylic acid 50% (paraffin)	14017601	1
Caelo	Salicylic acid 50% (paraffin)	14017601	2
Caelo	Salicylic acid 50% (paraffin)	1401761	1
Phoenix 2.07.14	Salicylic acid 50% (paraffin)	14107901	1
Fagron	Salicylic acid 50% (paraffin)	14108301	1
Caeo AHDE 25.11.2014	Salicylic acid 50% (paraffin)	14205101	1
Caelo	Salicylic acid 50% (paraffin)	14337101	1
sanacorp	Salicylic acid 50% (paraffin)	14337301	1
Caelo	Salicylic acid 50% (paraffin)	14337701	3
Caelo	Salicylic acid 50% (paraffin)	14337601	2
Caelo	Salicylic acid 50% (paraffin)	15231701	6
Noweda 26.08.15	Salicylic acid 50% (paraffin)	15231801	1
Caelo	Salicylic acid 50% (paraffin)	15247201	7
Caelo	Salicylic acid 50% (paraffin)	15277001	5
Caesar&Loretz GmbH	Salicylic acid 50% (paraffin)	15305901	1
Caelo	Salicylic acid 50% (paraffin)	15373001	6
Lamotte	Salicylic acid 50% (paraffin)	5375300	1
Henry Lamotte/Anzag	Salicylic acid 50% (paraffin)	5448100	1
Caelo	Salicylic acid 50% (paraffin)	17172701	1
Linden-Apo	Salicylic acid 50% (paraffin)	2615M-01046	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	5448100	1
AHD/ Henry Lamotte	Salicylic acid 50% (paraffin)	5763800	2
Caelo	Salicylic acid 50% (paraffin)	17105501	4
AHDE08.05.2015	Salicylic acid 50% (paraffin)	5763800	1
Lamotte	Salicylic acid 50% (paraffin)	5763800	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	5763800	2
Fiebig	Salicylic acid 50% (paraffin)	5763800	1
HenryLamotte	Salicylic acid 50% (paraffin)	5763800	1
Henry Lamotte Olis GmbH...	Salicylic acid 50% (paraffin)	5763800	1
Henry Lamotte/Sanacorp	Salicylic acid 50% (paraffin)	5809400	1
Henry Lamotte/Kehr	Salicylic acid 50% (paraffin)	5809400	1
noweda / lamotte	Salicylic acid 50% (paraffin)	5809400	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	5809400	1
Lamotte	Salicylic acid 50% (paraffin)	5809400	1
Caelo	Salicylic acid 50% (paraffin)	5809400	1
Lamotte/Anzag	Salicylic acid 50% (paraffin)	5809400	1
Sanacorp	Salicylic acid 50% (paraffin)	5809400	1
lamotte	Salicylic acid 50% (paraffin)	5809400	1
Henry Lamotte	Salicylic acid 50% (paraffin)	6001023001	1
Lamotte/PX	Salicylic acid 50% (paraffin)	5809400	1
Henry Lamotte	Salicylic acid 50% (paraffin)	5809400	1
Henry Lamotte / Sana	Salicylic acid 50% (paraffin)	5809400	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	6001023001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
HenryLamotte	Salicylic acid 50% (paraffin)	6001023001	1
Henry Lamotte/ Noweda	Salicylic acid 50% (paraffin)	6001023001	1
Lamotte/K	Salicylic acid 50% (paraffin)	6001023001	1
AHD/ Henry Lamotte	Salicylic acid 50% (paraffin)	6001023001	4
Lamotte	Salicylic acid 50% (paraffin)	6001023001	1
Henry Lamotte Olis GmbH...	Salicylic acid 50% (paraffin)	6001023001	1
Caelo	Salicylic acid 50% (paraffin)	6001023002	1
Lamotte	Salicylic acid 50% (paraffin)	6001023002	3
Henry Lamotte	Salicylic acid 50% (paraffin)	6001023002	1
Noweda	Salicylic acid 50% (paraffin)	6001023002	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	6001862002	1
Caelo	Salicylic acid 50% (paraffin)	AKD9	1
GEHE	Salicylic acid 50% (paraffin)	3A050318	1
Apomix / 11	Salicylic acid 50% (paraffin)	8A240616	1
GEHE	Salicylic acid 50% (paraffin)	4A050318	1
Caelo	Salicylic acid 50% (paraffin)	107105801	1
GEHE	Salicylic acid 50% (paraffin)	10A110817	1
GEHE	Salicylic acid 50% (paraffin)	13A041017	1
Euro OTC	Salicylic acid 50% (paraffin)	15305801	1
Caelo	Salicylic acid 50% (paraffin)	15305601	1
Phoenix	Salicylic acid 50% (paraffin)	15373101	1
Caelo	Salicylic acid 50% (paraffin)	15373201	1
Fagron	Salicylic acid 50% (paraffin)	15373301	1
Caelo	Salicylic acid 50% (paraffin)	16121203	1
Caelo	Salicylic acid 50% (paraffin)	16121301	4
Caelo	Salicylic acid 50% (paraffin)	16121404	1
Caelo	Salicylic acid 50% (paraffin)	1612301	1
Bombastus	Salicylic acid 50% (paraffin)	16197401	1
Sanacorp, 22.09.2016	Salicylic acid 50% (paraffin)	16197401	1
WE:Sanacorp 29.10.2016 ...	Salicylic acid 50% (paraffin)	1612401	1
Caeanol	Salicylic acid 50% (paraffin)	16197401	1
Noweda	Salicylic acid 50% (paraffin)	16276101	1
apomix/Noweda	Salicylic acid 50% (paraffin)	162796101	1
Caelo	Salicylic acid 50% (paraffin)	17017401	6
Allianz Healthcare	Salicylic acid 50% (paraffin)	17017401	1
Caelo	Salicylic acid 50% (paraffin)	171017501	1
Caelo	Salicylic acid 50% (paraffin)	17017501	2
Caelo	Salicylic acid 50% (paraffin)	170801SV/17172701	1
Caelo	Salicylic acid 50% (paraffin)	1710501	1
Ebert-Jacobi	Salicylic acid 50% (paraffin)	17105501	1
Caelo	Salicylic acid 50% (paraffin)	17105602	2
Caelo	Salicylic acid 50% (paraffin)	17105801	8
Gehe	Salicylic acid 50% (paraffin)	17105801	2
Caelo	Salicylic acid 50% (paraffin)	17105702	1
Caesar&Loretz GmbH	Salicylic acid 50% (paraffin)	17105801	1
Gehe	Salicylic acid 50% (paraffin)	17145401	1
Caelo	Salicylic acid 50% (paraffin)	17145401	5
Cealo	Salicylic acid 50% (paraffin)	17172701	1
Caelo	Salicylic acid 50% (paraffin)	17172702	2
Caelo	Salicylic acid 50% (paraffin)	17261301	1
SanacorpWE:18.10.17 EK...	Salicylic acid 50% (paraffin)	17172702	1
Caelo	Salicylic acid 50% (paraffin)	17261401	4
Caleo	Salicylic acid 50% (paraffin)	17261401	1
Caesar & Loretz/ Noweda	Salicylic acid 50% (paraffin)	17326201	1
Caelo	Salicylic acid 50% (paraffin)	18034301	3
GEHE	Salicylic acid 50% (paraffin)	18034301	1
Caelo	Salicylic acid 50% (paraffin)	6001023001	1
Lamotte	Salicylic acid 50% (paraffin)	6001862001	1
Henry Lamotte Oils GMBH	Salicylic acid 50% (paraffin)	6001862001/170307SV	1
Henry Lamotte Oils GMBH	Salicylic acid 50% (paraffin)	6001862001/161014SV	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte / Kehr	Salicylic acid 50% (paraffin)	6001023002	2
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	6001023002	1
Henry Lamotte Olis GmbH...	Salicylic acid 50% (paraffin)	6001862002	1
Henry Lamotte	Salicylic acid 50% (paraffin)	6001862002	1
Henry Lamotte/Noweda	Salicylic acid 50% (paraffin)	6001862002	1
Henry Lamotte	Salicylic acid 50% (paraffin)	6004854001	2
Henry Lamotte/Noweda	Salicylic acid 50% (paraffin)	6004854001	2
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	6001862002	1
henry Lamotte / Noweda	Salicylic acid 50% (paraffin)	6004854001	1
Henry Lamotte/Kehr	Salicylic acid 50% (paraffin)	6004854001	1
Apomix / 11	Salicylic acid 50% (paraffin)	4A200613	1
Fagron	Salicylic acid 50% (paraffin)	9D28	1
Henry Lamotte/Noweda	Salicylic acid 50% (paraffin)	6007886001	1
Phönix Gotha	Salicylic acid 50% (paraffin)	10333901	1
Apomix	Salicylic acid 50% (paraffin)	1A240113	1
Lamotte/Phönix	Salicylic acid 50% (paraffin)	6009168001	1
Caelo	Salicylic acid 50% (paraffin)	613I-01046	1
henry Lamotte / Noweda	Salicylic acid 50% (paraffin)	6007886001	1
Caelo	Salicylic acid 50% (paraffin)	11226801	4
Caelo	Salicylic acid 50% (paraffin)	11226601	1
Phönix	Salicylic acid 50% (paraffin)	11274301	2
Caelo	Salicylic acid 50% (paraffin)	11274301	1
Caelo	Salicylic acid 50% (paraffin)	11333602	1
Noweda	Salicylic acid 50% (paraffin)	11333701	1
Caelo	Salicylic acid 50% (paraffin)	11333801	7
Euro OTC	Salicylic acid 50% (paraffin)	11333801	1
Noweda	Salicylic acid 50% (paraffin)	11333801	1
Caelo	Salicylic acid 50% (paraffin)	12061201	9
Caelo	Salicylic acid 50% (paraffin)	12061301	1
Phönix, 30.10.2012	Salicylic acid 50% (paraffin)	12061301	1
Kehr	Salicylic acid 50% (paraffin)	12201901	1
Noweda	Salicylic acid 50% (paraffin)	12201901	1
Anzag 15.10.12	Salicylic acid 50% (paraffin)	12201901	1
Caelo	Salicylic acid 50% (paraffin)	12228301	13
Gehe	Salicylic acid 50% (paraffin)	12228301	1
Bombastus	Salicylic acid 50% (paraffin)	12228301	1
Caelo	Salicylic acid 50% (paraffin)	12305901	2
Caelo	Salicylic acid 50% (paraffin)	12228402	1
Gehe	Salicylic acid 50% (paraffin)	12305902	1
Caelo	Salicylic acid 50% (paraffin)	12305902	2
Sanacorp	Salicylic acid 50% (paraffin)	1305202	1
Caelo	Salicylic acid 50% (paraffin)	126061201	1
Caelo	Salicylic acid 50% (paraffin)	13052401	4
Phönix Gotha	Salicylic acid 50% (paraffin)	13052401	1
Caelo	Salicylic acid 50% (paraffin)	13052502	11
Caelo	Salicylic acid 50% (paraffin)	1309211046	1
Caelo	Salicylic acid 50% (paraffin)	1506211046	1
Caelo	Salicylic acid 50% (paraffin)	15373101	1
Caelo	Salicylic acid 50% (paraffin)	15231501	1
Caelo	Salicylic acid 50% (paraffin)	20310812-6	1
Caelo	Salicylic acid 50% (paraffin)	2002311046	1
Gehe	Salicylic acid 50% (paraffin)	4313M-07018	1
Defektur	Salicylic acid 50% (paraffin)	46/13	1
DAC 2005	Salicylic acid 50% (paraffin)	4960200	1
Henry Lamotte	Salicylic acid 50% (paraffin)	5253300	1
Henry Lamotte	Salicylic acid 50% (paraffin)	4960200	1
Henry Lamotte/Anzag	Salicylic acid 50% (paraffin)	5253300	1
Henry Lamotte/ Anzag	Salicylic acid 50% (paraffin)	5065800	1
Henry Lamotte Oils GmbH	Salicylic acid 50% (paraffin)	5253300	1
Phönix	Salicylic acid 50% (paraffin)	5375300	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Salicylic acid 50% (paraffin)	AKD2	1

- 90 666 spectra from 501 *Apo-Ident* customers from a total of 26 848 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Salicylic acid 50% (paraffin)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Salicylic acid 50% (paraffin)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	320	0	35 224
Type B	0	480	0	93 394
Type C	0	680	6	90 666

The substance/substance group *Salicylic acid 50% (paraffin)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.1250 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.7500 %)
Type C	100.0000 % (> 99.6586 %)	99.1254 % (> 98.6880 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30883	31903	1.99	146.79
30927	31727	6.42	156.27
31727	31727	0.00	154.63

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31901	31901	0.00	143.00
31902	31902	0.00	148.61
31903	31903	0.00	145.36
32176	31727	4.73	156.53
32520	32520	0.00	141.49

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Shea butter, refined
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30289-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Shea butter, refined

Special notes

When selecting the *Shea butter, refined* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Shea butter, refined*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Shea butter, ref. ..	13143819	31513	40	1412122
Caelo	Shea butter, ref. ..	14261107	31582	40	20140923
Caelo	Shea butter, ref. ..	15303004	31826	30	20151112
Caelo	Shea butter, ref. ..	15144903	32300	40	20150703
Caelo	Shea butter, ref. ..	15144903	32301	40	20150703

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 190 spectra of 5 reference samples from the substance/substance group *Shea butter, refined*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 354 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 150 spectra from 5 reference samples of the substance/substance group *Shea butter, refined*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Shea butter, refined	10133310	30289	40
Caelo	Shea butter, refined	13143819	31513 [†]	20
Caelo	Shea butter, refined	14261107	31582 [†]	20
Caelo	Shea butter, refined	15303004	31826 [†]	30
Caelo	Shea butter, refined	17085805	33470	40

- 93 724 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 90 spectra from 37 *Apo-Ident* customers from 61 batches from the substance/substance group *Shea butter, refined*.
- These include spectra of independent samples from 59 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Shea butter, refined	14033304	2
sanacorp	Shea butter, refined	14033310	1
Caelo	Shea butter, refined	14261102	1
Caelo	Shea butter, refined	14261107	2
Caelo	Shea butter, refined	14261108	4
Caelo	Shea butter, refined	15144902	1
Caelo	Shea butter, refined	16017203	3
Caelo	Shea butter, refined	17085806	1
Caelo	Shea butter, refined	10133305	2
Caelo	Shea butter, refined	13143807	1
Caelo	Shea butter, refined	1224210	1
Fagron	Shea butter, refined	13C04-N01	5
Caelo	Shea butter, refined	13143814	2
Caelo	Shea butter, refined	13143813	2
Caesar/Noweda	Shea butter, refined	13143813	2
Abis Pharma/Sanacorp	Shea butter, refined	10-500P/0714	1
Caelo	Shea butter, refined	14033312	2
Caelo	Shea butter, refined	14033306	1
Caelo	Shea butter, refined	14033313	1
Caelo	Shea butter, refined	14033314	1
Caelo	Shea butter, refined	14261110	1
Caelo	Shea butter, refined	14261111	1
Caelo	Shea butter, refined	15144903	1
Sheabutter Import Franz...	Shea butter, refined	14-3-0001	1
Caelo	Shea butter, refined	15144907	1
Sheabutter Import Franz...	Shea butter, refined	14314-18515	1
Caelo	Shea butter, refined	15303002	2
Caelo	Shea butter, refined	15303009	2
Caelo	Shea butter, refined	15303005	1
Caelo	Shea butter, refined	16017204	1
Caelo	Shea butter, refined	16017205	1
Caelo	Shea butter, refined	16017208	1
Caelo	Shea butter, refined	16017214	2
Caelo	Shea butter, refined	17085809	1
Caelo	Shea butter, refined	17031501	1
Caelo	Shea butter, refined	17085803	1
Sanacop WE:18.11.17 100...	Shea butter, refined	17085809	1
Caelo	Shea butter, refined	17287506	1
sanacorp WE:28.12.17 E...	Shea butter, refined	17085815	1
Euro OTC	Shea butter, refined	1802A0688	1
Farfalla Essen/ Bahnhof...	Shea butter, refined	FA0130140	1
DAC 2006	Shea butter, refined	11346509	1
Gehe	Shea butter, refined	11079802	1
Caelo	Shea butter, refined	11346509	1
Gehe	Shea butter, refined	11346509	1
Bahnhofapotheke Kempten...	Shea butter, refined	G16FGFI	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
ECOCERT/Sanacorp	Shea butter, refined	keine	2
Caelo	Shea butter, refined	11346514	1
Fa. Franz Thome	Shea butter, refined	1205355	2
Fagron	Shea butter, refined	121273	2
Fagron	Shea butter, refined	11J07-N08	1
Caelo	Shea butter, refined	12242404	1
Caelo	Shea butter, refined	12242406	1
Caelo	Shea butter, refined	12242410	2
Caelo	Shea butter, refined	12242412	1
Caelo	Shea butter, refined	12242416	2
Fagron	Shea butter, refined	12D16-N03	1
Fagron	Shea butter, refined	12G12-N14	1
Fagron	Shea butter, refined	12g12-n14	1
Caelo	Shea butter, refined	13143810	1
Fagron	Shea butter, refined	13c04-n01	2
Franz Thome	Shea butter, refined	400-15	1
Sheabutter-Import Franz...	Shea butter, refined	413-0001	1
Canea/Phönix	Shea butter, refined	SB1351M13	1
Wolfvital	Shea butter, refined	W68B058	1

- 91 262 spectra from 501 *Apo-Ident* customers from a total of 26 950 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Shea butter, refined* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Shea butter, refined* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	190	0	35 354
Type B	0	135	15	93 717
Type C	0	62	28	91 260

The substance/substance group *Shea butter, refined* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.8421 %)
Type B	100.0000 % (> 99.9762 %)	90.0000 % (> 88.0000 %)
Type C	100.0000 % (> 99.6680 %)	68.8889 % (> 65.5556 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to

exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31513	31513	0.00	28.37
31582	31513	7.58	31.81
31826	31513	12.29	33.42
32300	31513	10.70	31.42
32301	31513	4.39	28.26

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Sodium lactate solution 50%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30136-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Sodium lactate solution 50%; Natrium lacticum solutum 50%

Special notes

When selecting the *Sodium lactate solution 50%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Sodium lactate solution 50%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sodium lactate s...	12027310	30621	80	20120131
Caelo	Sodium lactate s...	14039503	31321	40	20140207
Caelo	Sodium lactate s...	16056201	32839	40	AR-17-FG-006366-01
Caelo	Sodium lactate s...	170977	33319	40	20170427*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 4 reference samples from the substance/substance group *Sodium lactate solution 50%*. These samples are listed above in the section *calibration samples*. The reference samples originate from 4 different batches.
- 35 344 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 173 spectra from 5 reference samples of the substance/substance group *Sodium lactate solution 50%*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Sodium lactate solution 50%	09258901	30136	20
Caelo	Sodium lactate solution 50%	11024003	30247	33
Caelo	Sodium lactate solution 50%	16056201	32586	40
Caelo	Sodium lactate solution 50%	16056201	32917	40
Caelo	Sodium lactate solution 50%	173153	34231	40

- 93 701 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 593 spectra from 250 *Apo-Ident* customers from 109 batches from the substance/substance group *Sodium lactate solution 50%*.
- These include spectra of independent samples from 105 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Sodium lactate solution 50%	13111704	11
Caelo	Sodium lactate solution 50%	13111705	5
Caelo	Sodium lactate solution 50%	13330302	6
Caelo	Sodium lactate solution 50%	14039511	4
Caelo	Sodium lactate solution 50%	14271905	1
Caelo	Sodium lactate solution 50%	14271901	17
Caesar & Loretz GmbH	Sodium lactate solution 50%	14271901	1
Caelo	Sodium lactate solution 50%	14353205	20
Caelo	Sodium lactate solution 50%	14353207	12
Caelo	Sodium lactate solution 50%	14353210	10
Caelo	Sodium lactate solution 50%	15135505	24
Caelo	Sodium lactate solution 50%	15251101	42
Caelo	Sodium lactate solution 50%	15375806	18
Caelo	Sodium lactate solution 50%	16056201	30
Caelo	Sodium lactate solution 50%	16321001	5
Caelo	Sodium lactate solution 50%	11238706	2
Caelo	Sodium lactate solution 50%	16321008	12
Caelo	Sodium lactate solution 50%	12294603	7
Caelo	Sodium lactate solution 50%	12294605	8
Caelo	Sodium lactate solution 50%	12294610	5
Caelo	Sodium lactate solution 50%	13267001	23
Caelo	Sodium lactate solution 50%	13267009	11
Caelo	Sodium lactate solution 50%	13330303	2
Caelo	Sodium lactate solution 50%	3041410	1
Caesar & Lorenz GmbH	Sodium lactate solution 50%	71014/1	1
Caelo	Sodium lactate solution 50%	10/AKD7	1
Gehe	Sodium lactate solution 50%	10094060	1
Caelo	Sodium lactate solution 50%	11024009	5
Caelo	Sodium lactate solution 50%	11/AKD16	1
Sanacorp	Sodium lactate solution 50%	11/AKD13	1
Caelo	Sodium lactate solution 50%	12294604	6
Fagron	Sodium lactate solution 50%	1311002383	1
Caelo	Sodium lactate solution 50%	13267002	2
Caelo	Sodium lactate solution 50%	13267003	9
Caelo	Sodium lactate solution 50%	13330305	11
Caelo	Sodium lactate solution 50%	13330307	9
Caelo	Sodium lactate solution 50%	13330306	1
Caelo	Sodium lactate solution 50%	13330308	1
Caelo	Sodium lactate solution 50%	133305	1
Fagron	Sodium lactate solution 50%	14/13	1
Caelo	Sodium lactate solution 50%	14039503	14
Euro OTC	Sodium lactate solution 50%		1
Euro OTC	Sodium lactate solution 50%	14039503	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium lactate solution 50%	14039507	3
Caelo	Sodium lactate solution 50%	14039508	2
Caelo	Sodium lactate solution 50%	1414Q03157	1
Caelo	Sodium lactate solution 50%	14182902	10
Caelo	Sodium lactate solution 50%	14182903	2
Caelo	Sodium lactate solution 50%	14182904	8
phoenix	Sodium lactate solution 50%	14271901	1
Phönix	Sodium lactate solution 50%	14271901	1
Fagron	Sodium lactate solution 50%	14271901	1
Caelo	Sodium lactate solution 50%	14271902	3
Caelo	Sodium lactate solution 50%	14271903	1
Caelo	Sodium lactate solution 50%	14271904	9
Caelo	Sodium lactate solution 50%	14271906	2
Caelo	Sodium lactate solution 50%	14271907	6
Euro OTC	Sodium lactate solution 50%	14271907	1
Gehe 12.1.15	Sodium lactate solution 50%	14271907	1
Fagron	Sodium lactate solution 50%	14271907	1
Caelo	Sodium lactate solution 50%	14271908	2
Caelo	Sodium lactate solution 50%	14353201	4
AHD 19.02.15	Sodium lactate solution 50%	14353202	1
Caelo	Sodium lactate solution 50%	14353202	8
Caelo	Sodium lactate solution 50%	14353204	2
Alliance health	Sodium lactate solution 50%	14353205	1
Noweda 18.05.2015	Sodium lactate solution 50%	14353205	1
Caesar & Loretz GmbH	Sodium lactate solution 50%	14353205	1
Caelo	Sodium lactate solution 50%	14353208	1
Phoenix	Sodium lactate solution 50%	14353210	1
Fagron	Sodium lactate solution 50%	1453202	1
Fagron	Sodium lactate solution 50%	14A21-810-295977	1
Fagron	Sodium lactate solution 50%	14A21-B10	2
Fagron	Sodium lactate solution 50%	14a21-b10-295977	1
Fagron	Sodium lactate solution 50%	14A21-B10-290492	2
Caelo	Sodium lactate solution 50%	14A21-B10-295977	1
Caelo	Sodium lactate solution 50%	15135501	1
Caelo	Sodium lactate solution 50%	15135508	2
Caelo	Sodium lactate solution 50%	1525110	1
Phoenix	Sodium lactate solution 50%	15251101	1
Caesar&Loretz	Sodium lactate solution 50%	15251101	1
Caesar & Loretz/Sanacorp	Sodium lactate solution 50%	15251101	1
Caelo	Sodium lactate solution 50%	15251102	10
Caelo	Sodium lactate solution 50%	15375805	1
Caesar & Loretz/Noweda	Sodium lactate solution 50%	15375806	1
Caesar&Loretz GmbH	Sodium lactate solution 50%	15375806	1
Caelo	Sodium lactate solution 50%	15375808	10
Ebert 02.08.2016	Sodium lactate solution 50%	15375808	1
Caelo	Sodium lactate solution 50%	1537806	1
Phoenix	Sodium lactate solution 50%	16056201	1
Caelo	Sodium lactate solution 50%	16056206	17
Sanacorp	Sodium lactate solution 50%	16056206	1
Caelo	Sodium lactate solution 50%	16056207	3
Henry Lamotte Oils GmbH	Sodium lactate solution 50%	16321001	1
Caelo	Sodium lactate solution 50%	16321002	6
Caelo	Sodium lactate solution 50%	16321003	3
Caelo	Sodium lactate solution 50%	16321004	1
Caelo	Sodium lactate solution 50%	16321009	2
Caelo	Sodium lactate solution 50%	17097704	14
Caelo	Sodium lactate solution 50%	17097701	2
Caelo	Sodium lactate solution 50%	17097708	1
Caelo	Sodium lactate solution 50%	17097709	2
Gehe	Sodium lactate solution 50%	17097709	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium lactate solution 50%	17315301	1
Caelo	Sodium lactate solution 50%	2014U-03157	1
Caelo	Sodium lactate solution 50%	11024001	1
Caelo	Sodium lactate solution 50%	507213157	1
Caelo	Sodium lactate solution 50%	11238701	4
Caelo	Sodium lactate solution 50%	11063006	1
Caelo	Sodium lactate solution 50%	11238707	2
Caelo	Sodium lactate solution 50%	11238712	2
Gehe	Sodium lactate solution 50%	11238712	1
Caelo	Sodium lactate solution 50%	12027301	1
Caelo	Sodium lactate solution 50%	12027303	3
Caelo	Sodium lactate solution 50%	12027310	3
phönix, 30.08.2012	Sodium lactate solution 50%	12027312	1
Caelo	Sodium lactate solution 50%	12180404	5
Caelo	Sodium lactate solution 50%	12180406	2
phönix	Sodium lactate solution 50%	12180404/8	1
Beiersdorf AG / Noweda	Sodium lactate solution 50%	12180406	1
Sanacorp	Sodium lactate solution 50%	12180407	1
Caelo	Sodium lactate solution 50%	12180407	5
Phönix	Sodium lactate solution 50%	12180407	1
Caelo	Sodium lactate solution 50%	1218407	1
Caelo	Sodium lactate solution 50%	12180408	2
Gehe	Sodium lactate solution 50%	12294603	1
Anzag 27.02.2013	Sodium lactate solution 50%	12294603	1
Kehr	Sodium lactate solution 50%	12294604	1
Noweda	Sodium lactate solution 50%	12294604	2
Caelo	Sodium lactate solution 50%	12294607	6
Caesar/Lorenz	Sodium lactate solution 50%	12294610	1
Fagron	Sodium lactate solution 50%	12J18-M04	5
Bergchemie	Sodium lactate solution 50%	13111704	1
Bombastus	Sodium lactate solution 50%	13111704	1
Phönix, 20.09.2013	Sodium lactate solution 50%	13111705	1
Noweda	Sodium lactate solution 50%	13267001	1
Caelo	Sodium lactate solution 50%	13267006	2
Caelo	Sodium lactate solution 50%	13267005	1
Gehe	Sodium lactate solution 50%	13267009	1
gehe	Sodium lactate solution 50%	13330302	1
Caelo	Sodium lactate solution 50%	13330304	2
Sanacorp	Sodium lactate solution 50%	14039508	1
Caelo	Sodium lactate solution 50%	14011301	1
Caesar&lorenz,Phoenix	Sodium lactate solution 50%	14182902	1
Caelo	Sodium lactate solution 50%	20310511-1	1
Caelo	Sodium lactate solution 50%	3111Q-03157	1
Caelo	Sodium lactate solution 50%		1
phönix	Sodium lactate solution 50%	9/AKD56	1
Phönix Gotha	Sodium lactate solution 50%	K40875122	1
Caelo	Sodium lactate solution 50%	N-260313-3	1

- 90 759 spectra from 501 *Apo-Ident* customers from a total of 26 901 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Sodium lactate solution 50%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Sodium lactate solution 50%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	172	1	93 701
Type C	0	579	14	90 759

The substance/substance group *Sodium lactate solution 50%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	99.4220 % (> 97.6879 %)
Type C	100.0000 % (> 99.6586 %)	97.6391 % (> 97.1332 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30621	32839	2.97	50.98
31321	32839	2.17	50.37
32839	32839	0.00	50.70
33319	33319	0.00	52.59

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Sodium laureth sulfate solution 27%
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30936-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Sodium laureth sulfate solution 27%; Natrii laurylethersulfatis solutio 27 %; Sodium laureth sulphate solution 27%; Solutio natrii laurylaethersulfurici 27 %

Special notes

When selecting the *Sodium laureth sulfate solution 27%* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Sodium laureth sulfate solution 27%*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sodium laureth s...	12315211	30936	50	20121029
Caelo	Sodium laureth s...	13069915	31114	39	20130402
Caelo	Sodium laureth s...	14049902	31349	40	8930
Caelo	Sodium laureth s...	161968	32696	30	20161010*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 159 spectra of 4 reference samples from the substance/substance group *Sodium laureth sulfate solution 27%*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 385 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 270 spectra from 8 reference samples of the substance/substance group *Sodium laureth sulfate solution 27%*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Sodium laureth sulfate solut...	10372703	30345	40
Caelo	Sodium laureth sulfate solut...	14049902	31349 [†]	20
Caelo	Sodium laureth sulfate solut...	161968	32696 [†]	10
Caelo	Sodium laureth sulfate solut...	16196801	32818	40
Caelo	Sodium laureth sulfate solut...	163071	33090	40
Caelo	Sodium laureth sulfate solut...	171210	33335	40
Caelo	Sodium laureth sulfate solut...	172233	33698	40
Caelo	Sodium laureth sulfate solut...	172233	33883	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 604 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 48 spectra from 21 *Apo-Ident* customers from 31 batches from the substance/substance group *Sodium laureth sulfate solution 27%*.
- These include spectra of independent samples from 29 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Sodium laureth sulfate solut...	12341104	1
Caelo	Sodium laureth sulfate solut...	13241101	1
Caesar & Loretz GmbH Hi...	Sodium laureth sulfate solut...	14049913	1
Caelo	Sodium laureth sulfate solut...	14197915	4
Caelo	Sodium laureth sulfate solut...	14355314	5
Caelo	Sodium laureth sulfate solut...	14355308	1
Phönix	Sodium laureth sulfate solut...	14355314	1
Caelo	Sodium laureth sulfate solut...	15145010	1
Caelo	Sodium laureth sulfate solut...	15145017	2
Caelo	Sodium laureth sulfate solut...	17121004	2
Caelo	Sodium laureth sulfate solut...	16354801	1
Caelo	Sodium laureth sulfate solut...	17121002	1
Caelo	Sodium laureth sulfate solut...	14355304	2
Caelo	Sodium laureth sulfate solut...	16307104	1
Sanacorp	Sodium laureth sulfate solut...	14049902	1
Caelo	Sodium laureth sulfate solut...	14049902	1
Caelo	Sodium laureth sulfate solut...	13241004	1
Caelo	Sodium laureth sulfate solut...	14197920	1
Caelo	Sodium laureth sulfate solut...	12315211	2
Caelo	Sodium laureth sulfate solut...	14049904	1
Caelo	Sodium laureth sulfate solut...	11322611	1
Caelo	Sodium laureth sulfate solut...	14355315	1
Caelo	Sodium laureth sulfate solut...	14355310	1
Caelo	Sodium laureth sulfate solut...	15145004	1
Caelo	Sodium laureth sulfate solut...	12044308	4
Caelo	Sodium laureth sulfate solut...	12315203	1
Caelo	Sodium laureth sulfate solut...	13069911	1
Caelo	Sodium laureth sulfate solut...	13069905	1
Caelo	Sodium laureth sulfate solut...	12044307	1
Caesar & Loretz GmbH, 4...	Sodium laureth sulfate solut...	12315211	1
Caelo	Sodium laureth sulfate solut...	13069903	1
Caelo	Sodium laureth sulfate solut...	13241001	1
Caelo	Sodium laureth sulfate solut...	14197916	1
Caelo	Sodium laureth sulfate solut...	2613M-01907	1

- 91 304 spectra from 501 *Apo-Ident* customers from a total of 26 979 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Sodium laureth sulfate solution 27%* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Sodium laureth sulfate solution 27%* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	159	0	35 385
Type B	0	253	17	93 604
Type C	0	48	0	91 304

The substance/substance group *Sodium laureth sulfate solution 27%* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2264 %)
Type B	100.0000 % (> 99.9761 %)	93.7037 % (> 92.5926 %)
Type C	100.0000 % (> 99.6587 %)	100.0000 % (> 87.5000 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30936	31349	4.08	29.71
31114	31349	4.63	27.92
31349	31349	0.00	25.35
32696	32696	0.00	22.73

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Soft ointment**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30896-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Soft ointment; Unguentum molle

Special notes

When selecting the *Soft ointment* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Soft ointment*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Soft ointment	12353301	30896	80	20121218
Henry Lamotte	Soft ointment	5215100	30924	50	20111111
Caelo	Soft ointment	14080501	31452	40	AR-15-FG-000142-01
Caelo	Soft ointment	153768	31984	40	20151103*
Caelo	Soft ointment	160833	32097	40	20160322*
Caelo	Soft ointment	163134	32871	40	20161116*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 290 spectra of 6 reference samples from the substance/substance group *Soft ointment*. These samples are listed above in the section *calibration samples*. The reference samples originate from 6 different batches.
- 35 254 spectra from a total of 834 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 1070 spectra from 25 reference samples of the substance/substance group *Soft ointment*.
- These include spectra of independent samples from 19 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Soft ointment	10377401	30291	40
Caelo	Soft ointment	12353301	30882	40
Caelo	Soft ointment	153768	31983	60
Caelo	Soft ointment	153768	31984 [†]	20
Caelo	Soft ointment	153770	31985	60
Caelo	Soft ointment	153770	31986	60
Caelo	Soft ointment	153772	31987	60
Caelo	Soft ointment	153772	31988	60
Caelo	Soft ointment	160830	32094	40
Caelo	Soft ointment	160831	32095	30
Caelo	Soft ointment	160832	32096	40
Caelo	Soft ointment	161463	32406	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Soft ointment	161464	32407	40
Caelo	Soft ointment	162452	32496	40
Caelo	Soft ointment	162453	32497	40
Caelo	Soft ointment	162454	32498	40
Caelo	Soft ointment	162455	32499	40
Caelo	Soft ointment	171969	33560	40
Caelo	Soft ointment	171971	33561	40
Caelo	Soft ointment	171972	33562	40
Caelo	Soft ointment	172631	33748	40
Caelo	Soft ointment	172632	33749	40
Caelo	Soft ointment	172632	33750	40
Caelo	Soft ointment	180021	34012	40
Caelo	Soft ointment	182012	34258	40

- 92 804 spectra from a total of 2146 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1474 spectra from 336 *Apo-Ident* customers from 292 batches from the substance/substance group *Soft ointment*.
- These include spectra of independent samples from 288 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Soft ointment	1097697	2
Caelo	Soft ointment	14007301	22
Caelo	Soft ointment	14007501	11
Caelo	Soft ointment	14007601	29
Caelo	Soft ointment	14202001	14
Anzag;27.10.14;3,28EUR	Soft ointment	14202301	2
Caelo	Soft ointment	14202301	33
Caelo	Soft ointment	14202302	7
Caelo	Soft ointment	14375601	15
Anzag;21.03.15;3,53EUR	Soft ointment	14375601	1
Caelo	Soft ointment	14375602	24
Caelo	Soft ointment	14375801	7
Caelo	Soft ointment	14375901	8
	Soft ointment	15059301	1
Caelo	Soft ointment	15059301	19
Caelo	Soft ointment	15204401	31
Caelo	Soft ointment	15204501	22
Caelo	Soft ointment	15376901	9

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Soft ointment	15377101	44
Caelo	Soft ointment	16083001	29
Caelo	Soft ointment	16083002	20
Caelo	Soft ointment	16083101	21
Caelo	Soft ointment	16083301	6
Caelo	Soft ointment	16245301	13
Caelo	Soft ointment	18002302	1
Bombastus	Soft ointment	284057	1
Henry Lamotte Oils	Soft ointment	5600500	3
HENRY LAMOTTE OILS GmbH	Soft ointment	6001368015	2
Anzag	Soft ointment	269470	1
Caelo	Soft ointment	11282101	3
Caelo	Soft ointment	12353401	11
Caelo	Soft ointment	13053602	6
Caelo	Soft ointment	13192501	27
Anzag	Soft ointment	267553	1
Caelo	Soft ointment	13192502	3
Caelo	Soft ointment	13192701	6
Caelo	Soft ointment	13375501	13
Lamotte	Soft ointment	5375800	1
Bombastus	Soft ointment	278047	1
Caelo	Soft ointment	13375402	11
Caelo	Soft ointment	13192601	36
Bombastus	Soft ointment	275422	4
Caelo	Soft ointment	16313401	18
GEHE	Soft ointment	16313501	1
Bombastus	Soft ointment	17000370	2
Caelo	Soft ointment	17076301	3
Ichthyol	Soft ointment	17076301	2
Caelo	Soft ointment	17076401	11
Noweda	Soft ointment	17076401	1
Caelo	Soft ointment	17076501	11
Caelo	Soft ointment	17076601	7
AHD	Soft ointment	17076601	1
Caelo	Soft ointment	17076701	6
Caelo	Soft ointment	17076801	2
Caelo	Soft ointment	17076901	3
Caelo	Soft ointment	17196901	3
Phönix	Soft ointment	17197101	1
Caelo	Soft ointment	17197101	4
Caelo	Soft ointment	17197201	1
GEHE	Soft ointment	17197102	1
GEHE	Soft ointment	17263101	1
Caelo	Soft ointment	17263101	4
Caelo	Soft ointment	17263201	13
Fagron	Soft ointment	17A03-T01-039997	5
Fagron	Soft ointment	17H28-T03-047827	6
Caelo	Soft ointment	18002101	3
Caelo	Soft ointment	18002303	1
Caelo	Soft ointment	18002501	2
GEHE	Soft ointment	18002501	1
Caelo	Soft ointment	18140502	1
Caelo	Soft ointment	8012015C	1
Caelo	Soft ointment	5081405	1
Caelo	Soft ointment	5011504	1
Caelo	Soft ointment	20170815-02542	1
Caelo	Soft ointment	7051510	1
Sanacorp	Soft ointment	100/AKD10	1
Caelo	Soft ointment	1097763	2
Caelo	Soft ointment	100/AKD7	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Soft ointment	101/AKD13	1
Caelo	Soft ointment	101/AKD12	1
caesar&lorenz GmbH,Phönix	Soft ointment	12194502	2
Anzag 18.04.14 EK:3,87	Soft ointment	1324-T02-006929	1
Gehe/04.12.14	Soft ointment	12353301	1
Caelo	Soft ointment	1324-T01-006925	1
Caelo	Soft ointment	13053801	19
Caelo	Soft ointment	1314I-02542	1
Caelo	Soft ointment	1314I2542	1
Fagron	Soft ointment	1324-t03-006933	1
Gehe am 25.02.14	Soft ointment	13375402	1
Caelo	Soft ointment	13375503	2
Caelo	Soft ointment	13387402	21
Berkel AHK	Soft ointment	13387402	1
Caelo	Soft ointment	13E13-T12-003670	1
Phoenix	Soft ointment	13387403	1
Fagron	Soft ointment	13I24-T01-006925	2
Euro OTC	Soft ointment	13E13-T13-003675	1
AHD 17.11.14 EK 5,51	Soft ointment	13I24-T01-006925	1
Fagron	Soft ointment	13I24-T02-006930	1
Fagron	Soft ointment	13I24-T01-006925	1
Fagron	Soft ointment	13I24-T03-006933	1
Fagron	Soft ointment	13I24-T04-006936	1
PHOENIX 10.03.2014	Soft ointment	14007302	1
Caelo	Soft ointment	14007302	2
Caelo	Soft ointment	14007402	3
	Soft ointment	14007601	1
Noweda	Soft ointment	14007601	1
Berkel AHK	Soft ointment	14007601	1
Caelo	Soft ointment	1402201	1
Caelo	Soft ointment	14080101	1
Caelo	Soft ointment	14080201	14
Euro OTC	Soft ointment	14080201	1
Alliance Healthcare	Soft ointment	14080201	1
Caelo	Soft ointment	14080301	5
Sanacorp	Soft ointment	14080301	1
Caelo	Soft ointment	14080501	9
Caesar &Loretz GmbH	Soft ointment	14080501	1
Caesar & Loretz	Soft ointment	14080501	1
Caelo	Soft ointment	14201801	2
Caelo	Soft ointment	14201901	7
Fagron	Soft ointment	14201901	1
Berkel AHK	Soft ointment	14202001	1
Sanacorp 19.01.15 EK 5,39	Soft ointment	14202001	1
Caelo	Soft ointment	14202011	1
Caelo	Soft ointment	14202201	25
Berkel AHK	Soft ointment	14202201	2
Sana 18.12.14 EK 5,39	Soft ointment	14202201	1
Gehe	Soft ointment	14202301	2
sanacorp	Soft ointment	14202301	1
Fagron	Soft ointment	14202301	1
Caelo	Soft ointment	14205001	1
Caelo	Soft ointment	1422001	1
Phönix	Soft ointment	14264301	1
Allianc 03.02.15	Soft ointment	14375601	1
Phönix	Soft ointment	14375601	1
Caesar&Loretz GmbH/Sana...	Soft ointment	14375602	1
AHDE 03.06.2015	Soft ointment	14375602	2
Allianz Healthcare 11.0...	Soft ointment	14375602	1
Euro OTC	Soft ointment	14375602	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Soft ointment	14375801	1
Caelo	Soft ointment	14375603	1
Caelo	Soft ointment	1437591	1
Fagron	Soft ointment	14E16-T03-012547	1
Fagron	Soft ointment	14e16-t03-012548	1
Fagron	Soft ointment	14E16-T03-012548	1
Cealo/Sanacorp	Soft ointment	15059301	1
Noweda	Soft ointment	15059301	1
Sanacorp	Soft ointment	15059301	1
Caelo	Soft ointment	15059401	9
Ichthyol-Gesellschaft C...	Soft ointment	15059401	1
Caelo	Soft ointment	15059402	2
Fagron	Soft ointment	15059601	1
Caelo	Soft ointment	15059601	3
Caelo	Soft ointment	15092014C	1
Caelo	Soft ointment	15059501	1
Caelo	Soft ointment	15204402	2
Caelo	Soft ointment	15204602	5
Caelo	Soft ointment	15204701	5
Caelo	Soft ointment	15204801	2
Caelo	Soft ointment	15376801	20
Sanacorp	Soft ointment	15376801	1
sanacoro	Soft ointment	15376801	1
Noweda	Soft ointment	15376901	2
Caelo	Soft ointment	15377001	19
Noweda	Soft ointment	15377001	2
InfectoPharm/Noweda	Soft ointment	15377001	1
Fagron	Soft ointment	15377001	1
Bombastus	Soft ointment	15377001	1
Frau Brinn-Holtz	Soft ointment	1577001	1
Fagron	Soft ointment	15A30	1
caesar&loretz GmbH,Phönix	Soft ointment	15959301	1
Fargron/ Noweda	Soft ointment	15A30-T02-018955	2
Fagron	Soft ointment	15A30-T02-018955	3
Caelo	Soft ointment	15e14-t01-022366	1
Fagron	Soft ointment	15E14-T02-022369	1
Fagron	Soft ointment	15F24-T16-023904	5
Fagron	Soft ointment	15F24T16023904	1
Fagron	Soft ointment	15f24-t16-023904	1
Fagron	Soft ointment	15I08-T01-0205888	1
Fagron	Soft ointment	15I08-T01-025888	1
Fagron	Soft ointment	15I08-T02-025891	1
Caelo	Soft ointment	16024401	1
Caelo	Soft ointment	15108t02025891	1
Caelo	Soft ointment	160620UM	1
Caelo	Soft ointment	1603001	1
phönix	Soft ointment	16083001	1
Phoenix	Soft ointment	16083002	1
phönix	Soft ointment	16083002	1
Caelo	Soft ointment	16083003	1
Wolf-Apotheke	Soft ointment	16083101	1
Bombastus	Soft ointment	16083101	1
Euro OTC	Soft ointment	16083101	1
Ebert & Jacobi 19.07...	Soft ointment	16083101	1
Caelo	Soft ointment	161111UM/16146103	1
Alte Apotheke	Soft ointment	16146103	1
Caelo	Soft ointment	16146103	3
Caelo	Soft ointment	1614620	1
Caelo	Soft ointment	16146201	5
Caelo	Soft ointment	16146301	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Soft ointment	16146301	1
Caelo	Soft ointment	16146302	1
Caelo	Soft ointment	16146401	7
Caelo	Soft ointment	16146502	12
Caesar & Loretz	Soft ointment	16245201	1
Caelo	Soft ointment	16146503	2
Caelo	Soft ointment	16245201	19
Phönix	Soft ointment	16245201	1
SanacorpWE:28.01.2017 ...	Soft ointment	16245201	1
Caelo	Soft ointment	16245202	4
Caelo	Soft ointment	16245401	10
Fagron	Soft ointment	16245401	1
Caelo	Soft ointment	16245402	5
Euro OTC	Soft ointment	16313401	1
Caelo	Soft ointment	16313402	1
Caelo	Soft ointment	16313501	17
Euro OTC	Soft ointment	16313501	3
Sanacorp	Soft ointment	16313501	1
Fagron	Soft ointment	16A21-T05-029597	5
Caelo	Soft ointment	16A21T05029597	1
Fagron	Soft ointment	16a21-t05-029597	1
Fagron	Soft ointment	16G18-T02-035247	3
Caelo	Soft ointment	16F03-T02-033989	1
Caelo	Soft ointment	16A21-T05-029597	1
Fagron	Soft ointment	16A21-T06-029597	1
Fagron	Soft ointment	16G18-T02-035248	3
phönix	Soft ointment	17076401	1
Caelo	Soft ointment	170509UM/16245202	1
Fagron	Soft ointment	17076401	1
Caeolo/Phönix	Soft ointment	17076401	1
Caelo	Soft ointment	17076402	1
Gehe	Soft ointment	17076402	1
Caelo	Soft ointment	170926UM/17076301	1
Fagron	Soft ointment	17A03-T01039997	1
Bombastus	Soft ointment	18000135	1
Caelo	Soft ointment	180320UM/17263302	1
Caelo	Soft ointment	20011702	1
Caelo	Soft ointment	23061403	1
Bombastus	Soft ointment	281897	1
Phönix Gotha	Soft ointment	286010	1
Bombastus	Soft ointment	284056	1
Fagron	Soft ointment	2514M-02542	1
Bombastus	Soft ointment	286012	9
Caelo	Soft ointment	286011	1
Bombastus	Soft ointment	288013	4
Bombastus	Soft ointment	288702	1
Bombastus	Soft ointment	288704	6
Phönix Gotha	Soft ointment	291331	1
Bombastus	Soft ointment	291333	3
BASF AG	Soft ointment	291333	1
Bombastus	Soft ointment	293336	3
Bomastus	Soft ointment	291333	1
Bombastus	Soft ointment	295313	3
Bombastus	Soft ointment	297077	9
Bombastus	Soft ointment	298069	1
Fagron	Soft ointment	299868	1
Bombastus	Soft ointment	298070	1
Bombastus	Soft ointment	299869	4
Bombastus	Soft ointment	300642	1
Caelo	Soft ointment	299869	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Soft ointment	302326	1
Bombastus	Soft ointment	302327	4
Bombastus	Soft ointment	304548	2
Bombastus	Soft ointment	305449	1
Bombastus	Soft ointment	305450	1
Bombastus	Soft ointment	305451	1
Bombastus	Soft ointment	309727	1
Caelo	Soft ointment	31031503	1
Henry Lamotte Oils GmbH...	Soft ointment	5448300	1
Caelo	Soft ointment	51123/1Sie	1
Bombastus	Soft ointment	3099727	1
Henry Lamotte Oils GmbH	Soft ointment	5600500	2
Henry Lamotte oils GmbH...	Soft ointment	5448300	1
Caelo	Soft ointment	5448300	1
Henry Lamotte oils GmbH...	Soft ointment	5691000	1
Lamotte/Alliance	Soft ointment	5691000	1
Henry Lamotte Oils GmbH	Soft ointment	5789000	2
AHD 02.09.14 EK 5,03	Soft ointment	5691800	1
Henry Lamotte Oils GmbH	Soft ointment	5691800	1
ahd/lamotte	Soft ointment	5789000	1
Lamotte/Kehr	Soft ointment	5789000	1
Henry Lamotte Oils GmbH...	Soft ointment	5789000	3
Lamotte	Soft ointment	5789000	2
Fagron	Soft ointment	5869200	1
Lamotte/Kehr	Soft ointment	5869200	1
Henry Lamotte/Anzag	Soft ointment	6001368001	1
Lamotte/Kehr	Soft ointment	6001368001	1
Henry Lamotte Oils GmbH	Soft ointment	6001368015	1
Henry Lamotte Oils GmbH...	Soft ointment	6001368015	2
Lamotte/Kehr	Soft ointment	6001368015	2
Henry Lamotte/Anzag	Soft ointment	6001368015	1
Henry Lamotte	Soft ointment	6001368015	2
Henry Lamotte Oils GmbH	Soft ointment	8002320003	1
Bombastus	Soft ointment	6862381	1
Henry Lamotte Oils / Gehe	Soft ointment	8002320003	1
Henry Lamotte/Geilenkir...	Soft ointment	8002320003	1
Lamotte/Kehr	Soft ointment	8002320004	1
Sanacorp	Soft ointment	94/AKD1	1
Sanacorp	Soft ointment	97/AKD4	1
Sanacorp	Soft ointment	96/AKD3	1
Phönix Gotha	Soft ointment	261107	1
Bombastus	Soft ointment	261106	1
Caelo	Soft ointment		1
Bombastus	Soft ointment	262274	1
Phönix Gotha	Soft ointment	268872	1
Bombastus	Soft ointment	269469	2
Bombastus	Soft ointment	269468	1
Caelo	Soft ointment	Caelo	1
Bombastus	Soft ointment	271430	1
Caelo	Soft ointment	612212542	1
Phönix Gotha	Soft ointment	271425	1
Caelo	Soft ointment	7071202	1
Caelo	Soft ointment	8012014C	1
Caelo	Soft ointment	11132901	1
Caelo	Soft ointment	10111201	1
Caelo	Soft ointment	11133001	2
Noweda	Soft ointment	11133001	1
Caelo	Soft ointment	11256301	1
Caelo	Soft ointment	11133002	1
Phönix	Soft ointment	11256301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Soft ointment	11256302	1
Caelo	Soft ointment	11256701	5
Noweda	Soft ointment	11256701	1
Caelo	Soft ointment	11282201	7
Caelo	Soft ointment	11282102	2
Noweda	Soft ointment	11282201	1
Caelo	Soft ointment	118282301	1
Caelo	Soft ointment	12045104	1
Caelo	Soft ointment	11282301	2
Caelo	Soft ointment	12045401	1
Caesar&loretz,Phoenix	Soft ointment	12045401	1
Caelo	Soft ointment	12045501	3
Caelo	Soft ointment	12045701	5
Caelo	Soft ointment	12045601	1
Caesar&loretz,Phoenix	Soft ointment	12045501	1
Phoenix	Soft ointment	12045701	1
Hustadt Apotheke	Soft ointment	12045701	1
Caelo	Soft ointment	12045702	5
Caelo	Soft ointment	12091401	5
Caesar&loretz,Phoenix	Soft ointment	12091303	1
Caelo	Soft ointment	120502/6	1
Caelo	Soft ointment	12091303	3
Gehe	Soft ointment	12091401	1
Anzag ek:5,12EUR 19.09.12	Soft ointment	12091402	1
Noweda	Soft ointment	12091402	1
Caelo	Soft ointment	12091402	6
Caelo	Soft ointment	12091403	1
Caelo	Soft ointment	12194401	3
Caelo	Soft ointment	12194501	8
Gehe	Soft ointment	12194501	2
Sanacorp	Soft ointment	12194501	1
Caelo	Soft ointment	12194502	3
Caelo	Soft ointment	12194601	13
Noweda	Soft ointment	12194601	1
Phönix	Soft ointment	12194601	1
Caelo	Soft ointment	12194701	2
phönix	Soft ointment	12194701	1
Caelo	Soft ointment	12353201	16
Großhandel	Soft ointment	12353201	1
Noweda	Soft ointment	12353201	1
Anzag 31.01.13 ek: 5,10EUR	Soft ointment	12353201	1
Caelo	Soft ointment	12353301	10
Anzag 23.01.2013	Soft ointment	12353301	1
fargon/Sanacorp	Soft ointment	12353301	1
Caelo	Soft ointment	12353302	2
Bombastus	Soft ointment	12353301	1
Caelo	Soft ointment	1253201	1
Caelo	Soft ointment	1294502	1
Fagron	Soft ointment	12e23-t01	1
Fagron	Soft ointment	12E23-T01	1
Fagron	Soft ointment	12G09-T02	1
Caelo	Soft ointment	130221/1	1
Kehr	Soft ointment	13053601	1
Caelo	Soft ointment	13053601	17
Caesar & Loretz GmbH, 4...	Soft ointment	13053601	1
Caesar u Loretz	Soft ointment	13053601	1
Gehe	Soft ointment	13053601	2
Caelo	Soft ointment	13053701	1
Großhandel	Soft ointment	13053801	1
Phönix	Soft ointment	13053801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
20; 27.08.13	Soft ointment	13053801	1
Caelo	Soft ointment	13053901	6
Caelo	Soft ointment	1319201	1
Gehe 10.10.13	Soft ointment	13192601	1
Phoenix	Soft ointment	13192601	1
gehe 04.12.13	Soft ointment	13192601	1
Sanacorp	Soft ointment	13192701	1
Caelo	Soft ointment	13192801	2
Caelo	Soft ointment	13192901	18
Sanacorp	Soft ointment	13192901	1
Phönix	Soft ointment	13192901	1
Caelo	Soft ointment	13192902	3
Euro OTC	Soft ointment	13192902	1
Sanacorp	Soft ointment	13375402	1
Fagron	Soft ointment	13387402	1
Caelo	Soft ointment	13387403	4
Gehe am 04.02.2014	Soft ointment	13387402	1
Fagron	Soft ointment	13A08-T12-000186	1
Fagron	Soft ointment	13A08-T13-000189	1
Fagron	Soft ointment	13C04-T04-0002	1
Fagron	Soft ointment	13C04-T04-001709	3
Fagron	Soft ointment	13C04-T04-001710	3
Fagron	Soft ointment	13C04-T05	2
Fagron	Soft ointment	13E13-T12-003670	3
Fagron	Soft ointment	13E13-T12-003669	1
Fagron	Soft ointment	13E13-T14-003679	2
Phoenix	Soft ointment	14007601	1
Phönix	Soft ointment	14201801	1
Caelo	Soft ointment	1420220	1
Caelo	Soft ointment	18121206	1
Caelo	Soft ointment	1912A-02542	1
Caelo	Soft ointment	20230712-4	1
Caelo	Soft ointment	2513A-02542	1
Bombastus	Soft ointment	275424	1
Bombastus	Soft ointment	278048	2
Bombastus	Soft ointment	276854	1
Phönix Gotha	Soft ointment	278046	1
Bombastus	Soft ointment	280166	2
Alliance Healthcare/Hed...	Soft ointment	281897	2
Caelo	Soft ointment	284057	1
Caelo	Soft ointment	3213E2542	1
Fagron	Soft ointment	4113A-02542	1
Henry Lamotte,Anzag	Soft ointment	5189900	2
Henry Lamotte,Phönix	Soft ointment	5067500	1
H. Lamotte	Soft ointment	5067500	1
Lamotte	Soft ointment	5215100	1
HenryLamotte/Phoenix	Soft ointment	5189900	1
Henry Lamotte	Soft ointment	5067500	1
Finze	Soft ointment	5215100	1
Henry Lamotte/GEHE	Soft ointment	5215100	1
Henry Lam./ Sanacorp	Soft ointment	5215100	1
lamotte/Noweda	Soft ointment	5265300	2
Henry Lamotte Oils GmbH	Soft ointment	5317100	3
lamotte7anzag	Soft ointment	5294600	1
Henry Lamotte/GEHE	Soft ointment	5317100	1
Henry Lamotte Oils	Soft ointment	5317100	1
Fiebig 32.01.13 ek:4,84EUR	Soft ointment	5317100	1
H. Lamotte	Soft ointment	5317100	1
Henry Lamotte	Soft ointment	5317100	1
H.Lamotte / Anzag	Soft ointment	5317100	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Soft ointment	5317100	1
Henry Lamotte oils GmbH	Soft ointment	5321600	2
Henra Lamotte	Soft ointment	5321600	1
Gehe	Soft ointment	5375800	1
Lamotte/ Anzag	Soft ointment	5375800	1
Henry Lamotte/Anzag	Soft ointment	5448300	1
Lamotte/ Anzag	Soft ointment	5496100	1
Henry Lamotte Oils GmbH...	Soft ointment	5496100	2
Henry Lamotte Oils GMBH	Soft ointment	5496100	1
Henry Lamotte/GEHE	Soft ointment	5375800	1
Henry Lamotte Oils GmbH	Soft ointment	5496100	1
Caelo	Soft ointment	I3192502	1
Noweda	Soft ointment	60089901	1
Lamotte/Krieger/20.05.2...	Soft ointment	5496100	1
Henry Lamotte Oils/ Gehe	Soft ointment	61580	1
Caelo	Soft ointment	92810389	1

- 89 878 spectra from 501 *Apo-Ident* customers from a total of 26 721 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Soft ointment* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Soft ointment* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	290	0	35 254
Type B	0	1039	31	92 804
Type C	0	1411	63	89 878

The substance/substance group *Soft ointment* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 97.9310 %)
Type B	100.0000 % (> 99.9761 %)	97.1028 % (> 96.8224 %)
Type C	100.0000 % (> 99.6586 %)	95.7259 % (> 95.5224 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30896	31984	5.31	101.10
30924	31984	10.37	115.86
31452	31452	0.00	93.73
31984	31984	0.00	98.93
32097	32097	0.00	87.38
32871	32871	0.00	95.99

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Softisan® 378**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30887-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Softisan® 378; Adeps solidus 378

Special notes

When selecting the *Softisan® 378* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Softisan® 378*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Softisan® 378	12K27-N05	30887	60	20121204
Fagron	Softisan® 378	13D26-N12	31117	40	1407632
Caelo	Softisan® 378	13370104	31472	40	20131119
Caelo	Softisan® 378	14301114	32298	30	20141118
Fagron	Softisan® 378	15G20-B10-326361	32964	40	1704407

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 210 spectra of 5 reference samples from the substance/substance group *Softisan*[®] 378. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 334 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 170 spectra from 5 reference samples of the substance/substance group *Softisan*[®] 378.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	<i>Softisan</i> [®] 378	10F09-N12	30295	40
Caelo	<i>Softisan</i> [®] 378	14301114	32298 [†]	10
Caelo	<i>Softisan</i> [®] 378	14301114	32299	40
Caelo	<i>Softisan</i> [®] 378	162686	33078	40
Caelo	<i>Softisan</i> [®] 378	162686	33101	40

- 93 704 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 136 spectra from 46 *Apo-Ident* customers from 77 batches from the substance/substance group *Softisan*® 378.
- These include spectra of independent samples from 72 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
THC Pharma / Phoenix	Softisan® 378	103897	2
Bionorica/Fiebig	Softisan® 378	109779	1
Caelo	Softisan® 378	16010408	1
Fagron	Softisan® 378	15G20-B10	1
Sino Phyto	Softisan® 378	2012060444	1
Fagron	Softisan® 378	17B16-B06-331749	1
GEHE	Softisan® 378	17B16-B06	1
Caelo	Softisan® 378	18011302	1
Allianc. Health. / THC ...	Softisan® 378	92456	1
Aallianc. Healthc. / TH...	Softisan® 378	98561	7
THC Pharma/Phönix	Softisan® 378	95076	1
Bionorica/Kehr	Softisan® 378	103897	1
Bionorica	Softisan® 378	107075	1
Bionorica/Kehr	Softisan® 378	107075	1
Bionorica Ethics	Softisan® 378	109779	1
Bionorica	Softisan® 378	109779	1
Caelo	Softisan® 378	13370104	3
Bionorica SE	Softisan® 378	107075	1
Caelo	Softisan® 378	3/2017	1
Fagron	Softisan® 378	514E-02911	1
Caelo	Softisan® 378	14301104	2
Caelo	Softisan® 378	14301113	1
Fagron	Softisan® 378	14F11-B04	3
Fagron	Softisan® 378	14F11-B04-296083	3
Fagron	Softisan® 378	14F11-B04-299298	3
Fagron	Softisan® 378	15G20-B10-316518	3
Caelo	Softisan® 378	15G20B10320309	1
Fagron	Softisan® 378	15G20-B10-525272	1
Fagron	Softisan® 378	15G20-B10-525116	1
Fagron	Softisan® 378	15G20-B10-525379	1
Caelo	Softisan® 378	16268605	1
Caelo	Softisan® 378	16010402	1
Caelo	Softisan® 378	15gb10326361	1
Caelo	Softisan® 378	17247103	1
Fagron	Softisan® 378	17B16-B06-525441	1
THC Pharma	Softisan® 378	2013091014	2
Caelo	Softisan® 378	22062015B	1
Fagron	Softisan® 378	311662	2
Fagron	Softisan® 378	505717	1
Bionorica/Gehe	Softisan® 378	75850	3
Gehe Kassel 09.02.2015	Softisan® 378	75850	1
Caelo	Softisan® 378	1081408	1
Caelo	Softisan® 378	4041501	1
Fagron	Softisan® 378	12I24-N01	1
Caelo	Softisan® 378	13370112	1
Caelo	Softisan® 378	13370102	2
Fagron	Softisan® 378	13D26-N12	12
Fagron	Softisan® 378	13H22-N01	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Softisan [®] 378	13H22-N01	1
Fagron	Softisan [®] 378	14F11B04-296083	1
Fagron	Softisan [®] 378	15G20-B10-310231	1
Fagron	Softisan [®] 378	15G20-10-316518	1
Fagron	Softisan [®] 378	15G20-B10-311893	1
Caelo	Softisan [®] 378	16021507	1
Thc	Softisan [®] 378	2012060444	1
Alliance Health.14.08.2...	Softisan [®] 378		1
Sanacorp	Softisan [®] 378	3/AKD1	1
Caelo	Softisan [®] 378	28031401	1
Caelo	Softisan [®] 378	57657	1
Fagron	Softisan [®] 378	303624	2
Fagron	Softisan [®] 378	11520	2
Caelo	Softisan [®] 378	73059	1
Caelo	Softisan [®] 378	3071201	1
Caelo	Softisan [®] 378	10310003	1
Caelo	Softisan [®] 378	11121304	1
Caelo	Softisan [®] 378	10310008	1
Fagron	Softisan [®] 378	12/24-N01	1
Caelo	Softisan [®] 378	13012703	1
Gehe	Softisan [®] 378	12D02-N07	1
Caelo	Softisan [®] 378	13012707	5
Fagron	Softisan [®] 378	12K27-N05	1
Caelo	Softisan [®] 378	11359302	1
Sanacorp	Softisan [®] 378	12A18-N11	1
Caelo	Softisan [®] 378	13012704	1
Fagron	Softisan [®] 378	13B11-N01	1
Fagron	Softisan [®] 378	13D26-N02	1
Fagron	Softisan [®] 378	13H22-N01	1
Caelo	Softisan [®] 378	14111312	1
Caelo	Softisan [®] 378	14301103	2
Caelo	Softisan [®] 378	14301112	1
Caelo	Softisan [®] 378	15G20-B10-310231	1
THC Pharma	Softisan [®] 378		1
Fagron	Softisan [®] 378	209602	4
Noweda	Softisan [®] 378	2012060444	1
Caelo	Softisan [®] 378	21031310	1
Caelo	Softisan [®] 378	209602	1
Caelo	Softisan [®] 378	27091306	1
Fagron	Softisan [®] 378	011520	1

- 91 216 spectra from 501 *Apo-Ident* customers from a total of 26 932 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Softisan[®] 378* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Softisan[®] 378* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	210	0	35 334
Type B	0	167	3	93 618
Type C	0	108	28	91 210

The substance/substance group *Softisan*[®] 378 can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.1429 %)
Type B	100.0000 % (> 99.9763 %)	98.2353 % (> 96.4706 %)
Type C	100.0000 % (> 99.6595 %)	79.4118 % (> 77.2059 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30887	31117	7.70	14.59
31117	31117	0.00	13.20
31472	31117	3.80	12.84
32298	31117	0.66	13.30
32964	32964	0.00	16.18

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Sorbitan- and glycerol monooleate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31574-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Sorbitan- and glycerol monooleate

Special notes

When selecting the *Sorbitan- and glycerol monooleate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Sorbitan- and glycerol monooleate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sorbitan- and gl...	14202402	31574	40	AR-15-FG-009320-01
Caelo	Sorbitan- and gl...	15122301	31726	40	20150420

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Sorbitan- and glycerol monooleate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 212 spectra from 7 reference samples of the substance/substance group *Sorbitan- and glycerol monooleate*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Sorbitan- and glycerol monoo...	14202402	31574 [†]	20
Caelo	Sorbitan- and glycerol monoo...	15122301	31726 [†]	20
Caelo	Sorbitan- and glycerol monoo...	16234301	33454	40
Caelo	Sorbitan- and glycerol monoo...	16234301	33455	40
Caelo	Sorbitan- and glycerol monoo...	17132202	33936	40
Caelo	Sorbitan- and glycerol monoo...	182151	34252	40
Caelo	Sorbitan- and glycerol monoo...	16234301	not required	12

- 93 662 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2 spectra from 2 *Apo-Ident* customers from 1 batches from the substance/substance group *Sorbitan- and glycerol monooleate*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Apomix	Sorbitan- and glycerol monoo...	20515	2

- 91 350 spectra from 501 *Apo-Ident* customers from a total of 27 008 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Sorbitan- and glycerol monooleate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Sorbitan- and glycerol monooleate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	212	0	93 662
Type C	0	2	0	91 350

The substance/substance group *Sorbitan- and glycerol monooleate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.1698 %)
Type C	100.0000 % (> 99.6630 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31574	31574	0.00	30.84
31726	31574	6.80	33.23

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Sorbitan trioleate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31031-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Sorbitan trioleate; Sorbitani trioleas

Special notes

When selecting the *Sorbitan trioleate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Sorbitan trioleate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sorbitan trioleate	12171410	31031	39	1407273
Caelo	Sorbitan trioleate	13317202	31130	40	20131017
Caelo	Sorbitan trioleate	13317212	31459	40	20131017

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 119 spectra of 3 reference samples from the substance/substance group *Sorbitan trioleate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 425 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 157 spectra from 5 reference samples of the substance/substance group *Sorbitan trioleate*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Sorbitan trioleate	10109108	30316	17
Caelo	Sorbitan trioleate	13317212	31459 [†]	20
Caelo	Sorbitan trioleate	15337406	33441	40
Caelo	Sorbitan trioleate	17056001	33461	40
Caelo	Sorbitan trioleate	180779	34080	40

- 93 717 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 44 spectra from 20 *Apo-Ident* customers from 25 batches from the substance/substance group *Sorbitan trioleate*.
- These include spectra of independent samples from 24 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Sorbitan trioleate	11284607	1
Caelo	Sorbitan trioleate	12171408	3
Gehe	Sorbitan trioleate	12171408	1
Caesar/Lorenz	Sorbitan trioleate	13317210	1
Caelo	Sorbitan trioleate	13317212	6
Caelo	Sorbitan trioleate	13317213	1
Caelo	Sorbitan trioleate	14218202	1
Caelo	Sorbitan trioleate	15064901	1
Caelo	Sorbitan trioleate	14218205	1
Caelo	Sorbitan trioleate	15064903	2
Caesar&Loretz GmbH	Sorbitan trioleate	15064903	1
Caelo	Sorbitan trioleate	15064908	1
Caelo	Sorbitan trioleate	15064909	1
Caelo	Sorbitan trioleate	15182502	2
Caelo	Sorbitan trioleate	15337404	2
Caelo	Sorbitan trioleate	15337403	1
Caelo	Sorbitan trioleate	15337405	3
GEHE	Sorbitan trioleate	15337404	1
Caesar & Lorenz GmbH	Sorbitan trioleate	15337404	1
Caelo	Sorbitan trioleate	17056001	1
Caelo	Sorbitan trioleate	15337406	1
Caelo	Sorbitan trioleate	17056004	1
GEHE	Sorbitan trioleate	17056001	1
Caelo	Sorbitan trioleate	11284603	1
Caelo	Sorbitan trioleate	11284608	1
Caesar &Loretz GmbH	Sorbitan trioleate	12171405	1
Caelo	Sorbitan trioleate	12171405	1
Caelo	Sorbitan trioleate	1217408	1
Caelo	Sorbitan trioleate	11284601	1
Caelo	Sorbitan trioleate	13317204	2
Caelo	Sorbitan trioleate	13317209	1

- 91 308 spectra from 501 *Apo-Ident* customers from a total of 26 984 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Sorbitan trioleate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Sorbitan trioleate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The

following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	119	0	35 425
Type B	0	157	0	93 717
Type C	0	44	0	91 308

The substance/substance group *Sorbitan trioleate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 94.9580 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.1783 %)
Type C	100.0000 % (> 99.6588 %)	100.0000 % (> 86.3636 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31031	31031	0.00	26.68
31130	31031	7.53	30.78
31459	31031	5.38	29.61

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Spearmint oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30668-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Spearmint oil; Menthae crispae aetheroleum; Oleum menthae crispae

Special notes

When selecting the *Spearmint oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Spearmint oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Spearmint oil	11346805	30668	40	20111206
Caelo	Spearmint oil	13363904	31303	40	1407357
Caelo	Spearmint oil	13363904	31419	40	20131118

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Spearmint oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 424 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Spearmint oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Spearmint oil	13363904	31419 [†]	20
Caelo	Spearmint oil	161284	32205	40
Caelo	Spearmint oil	172396	33912	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 12 spectra from 9 *Apo-Ident* customers from 11 batches from the substance/substance group *Spearmint oil*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Spearmint oil	14147608	1
Noweda	Spearmint oil	11346805	1
Caelo	Spearmint oil	16128409	1
Caelo	Spearmint oil	12206913	1
Caelo	Spearmint oil	17239605	1
Caelo	Spearmint oil	16128404	1
Caelo	Spearmint oil	25062015B	1
Caelo	Spearmint oil	16128402	1
Caelo	Spearmint oil	12206910	1
Caelo	Spearmint oil	16128405	1
Caelo	Spearmint oil	13093402	1
Sanacorp	Spearmint oil	13093402	1

- 91 340 spectra from 501 *Apo-Ident* customers from a total of 26 998 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Spearmint oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Spearmint oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	100	0	93 774
Type C	0	12	0	91 340

The substance/substance group *Spearmint oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6593 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30668	31303	11.63	39.88
31303	31303	0.00	41.05
31419	31303	7.14	43.23

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Spirit of ether**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30707-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Spirit of ether; Spiritus aethereus

Special notes

When selecting the *Spirit of ether* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Spirit of ether*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Spirit of ether	12074603	30707	40	20120320
Caelo	Spirit of ether	13131805	31207	40	1403111

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *Spirit of ether*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Spirit of ether*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Spirit of ether	15442201	32566	40
Caelo	Spirit of ether	15442203	33218	40
Caelo	Spirit of ether	181156	34065	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 10 spectra from 4 *Apo-Ident* customers from 7 batches from the substance/substance group *Spirit of ether*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Spirit of ether	13131805	2
Gehe	Spirit of ether	15442203	1
Caelo	Spirit of ether	12074603	2
Caelo	Spirit of ether	13131807	1
Caelo	Spirit of ether	11375509	1
Caelo	Spirit of ether	14339702	2
Caelo	Spirit of ether	11375512	1

- 91 342 spectra from 501 *Apo-Ident* customers from a total of 27 002 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Spirit of ether* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Spirit of ether* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	120	0	93 754
Type C	0	9	1	91 342

The substance/substance group *Spirit of ether* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6594 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30707	31207	5.87	129.94
31207	31207	0.00	140.59

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Spirit soap**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30546-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Spirit soap; Spiritus saponatus

Special notes

When selecting the *Spirit soap* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Spirit soap*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Spirit soap	12247110	31153	40	1403401
Caelo	Spirit soap	13278409	31393	40	20130829
Caelo	Spirit soap	173045	33893	40	20171204*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Spirit soap*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 160 spectra from 4 reference samples of the substance/substance group *Spirit soap*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Spirit soap	11265306	30546	40
Caelo	Spirit soap	11265302	30588	40
Caelo	Spirit soap	16130209	33443	40
Caelo	Spirit soap	16130209	33456	40

- 93 714 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 97 spectra from 40 *Apo-Ident* customers from 42 batches from the substance/substance group *Spirit soap*.
- These include spectra of independent samples from 40 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Spirit soap	12247108	7
Caelo	Spirit soap	13278409	6
Caelo	Spirit soap	12247109	2
Caelo	Spirit soap	12247110	7
Caelo	Spirit soap	6081516	1
Caelo	Spirit soap	13278403	8
Anzag 8.10.14 EK 4.88	Spirit soap	13278403	1
Caelo	Spirit soap	12121402	1
Caelo	Spirit soap	13278402	4
Caelo	Spirit soap	13278405	3
Caelo	Spirit soap	9031501	1
Caelo	Spirit soap	13278408	2
Anzag	Spirit soap	13278409	1
Caelo	Spirit soap	14212703	1
Caelo	Spirit soap	14212705	1
Caelo	Spirit soap	14212706	6
Caelo	Spirit soap	14212710	1
Caelo	Spirit soap	14212714	1
Caelo	Spirit soap	14212716	1
Caelo	Spirit soap	14212712	1
Caelo	Spirit soap	14364902	3
Caelo	Spirit soap	14364904	3
Caelo	Spirit soap	14364901	1
Caesar&Lorenz	Spirit soap	14364906	1
Caelo	Spirit soap	14364906	4
Caelo	Spirit soap	14364907	1
Caelo	Spirit soap	14364909	1
Caelo	Spirit soap	16069902	1
Caelo	Spirit soap	16130201	3
Fagron	Spirit soap	16130201	1
Caelo	Spirit soap	16130202	2
Caelo	Spirit soap	16130204	3
Caelo	Spirit soap	16130209	1
GEHE	Spirit soap	17304502	1
Caelo	Spirit soap	17304504	1
Caelo	Spirit soap	10121303	1
Caelo	Spirit soap	17304502	1
Sanacorp	Spirit soap	12247102	1
Caelo	Spirit soap	3091306	1
Caelo	Spirit soap	9172002	1
Caelo	Spirit soap	12247102	1
Caelo	Spirit soap	12247104	2
Caelo	Spirit soap	12247105	1
Caelo	Spirit soap	12247107	1
Caelo	Spirit soap	14364905	1
Caelo	Spirit soap	2712Q-02323	1
Caelo	Spirit soap	27072013	1
Caelo	Spirit soap	3913U-02323	1

- 91 255 spectra from 501 *Apo-Ident* customers from a total of 26 968 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Spirit soap* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Spirit soap* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	160	0	93 714
Type C	0	95	2	91 255

The substance/substance group *Spirit soap* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 96.2500 %)
Type C	100.0000 % (> 99.6586 %)	97.9381 % (> 94.8454 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31153	31153	0.00	40.21
31393	31153	2.52	41.42
33893	33893	0.00	46.14

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Spruce needle oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31414-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Spruce needle oil; Oleum pini sibiricum; Piceae aetheroleum

Special notes

When selecting the *Spruce needle oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Spruce needle oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Spruce needle oil	13152004	31414	30	1412281
Caelo	Spruce needle oil	13152004	31618	40	20130617
Caelo	Spruce needle oil	14348607	31787	40	20141202

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 110 spectra of 3 reference samples from the substance/substance group *Spruce needle oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 434 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 130 spectra from 5 reference samples of the substance/substance group *Spruce needle oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Spruce needle oil	13152004	31414 [†]	10
Caelo	Spruce needle oil	13152004	31618 [†]	20
Caelo	Spruce needle oil	14348607	31787 [†]	20
Caelo	Spruce needle oil	163511	33095	40
Caelo	Spruce needle oil	171966	33682	40

- 93 744 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 7 spectra from 5 *Apo-Ident* customers from 6 batches from the substance/substance group *Spruce needle oil*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Spruce needle oil	13389508	1
Caelo	Spruce needle oil	15418210	1
Caelo	Spruce needle oil	14348607	1
Caelo	Spruce needle oil	17196602	2
Taoasis	Spruce needle oil	22336-122331	1
Nachprüfung	Spruce needle oil	2914A-02006	1

- 91 345 spectra from 501 *Apo-Ident* customers from a total of 27 003 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Spruce needle oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Spruce needle oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	110	0	35 434
Type B	0	130	0	93 744
Type C	0	5	2	91 345

The substance/substance group *Spruce needle oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 94.5455 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.3846 %)
Type C	100.0000 % (> 99.6598 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31414	31414	0.00	15.53
31618	31414	2.64	16.48
31787	31414	6.14	16.58

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Squalane**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31327-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Squalane; Squalanum

Special notes

When selecting the *Squalane* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Squalane*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Squalane	13202603	31327	40	1504359
Caelo	Squalane	14330704	32075	40	9465
Caelo	Squalane	16030504	32302	40	20160205
Caelo	Squalane	171380	33317	40	20170609*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Squalane*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Squalane*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Squalane	13202603	31328	40
Caelo	Squalane	13202603	31329	40
Caelo	Squalane	17138001	33422	40

- 93 754 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 6 spectra from 2 *Apo-Ident* customers from 5 batches from the substance/substance group *Squalane*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Squalane	16030504	1
Spinnrad/Noweda	Squalane	116001301	1
Euro OTC	Squalane	316001302	1
Testsubstanz Saficalcan	Squalane	3104	1
Kaden	Squalane		2

- 91 346 spectra from 501 *Apo-Ident* customers from a total of 27 005 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Squalane* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Squalane* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	120	0	93 754
Type C	0	3	3	91 346

The substance/substance group *Squalane* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6600 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31327	31327	0.00	12.12
32075	32075	0.00	13.17
32302	32075	2.82	11.53
33317	33317	0.00	13.44

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Syrups**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30753-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Syrups; Althae sirupus (konserviert); Cherry syrup (conserved); Ipecacuanha syrup (NRF 19.1); Marshmallow syrup (conserved); Raspberry syrup (conserved); Ribwort plantain syrup (conserved); Sirupus althaeae (konserviert); Sirupus cerasi (konserviert); Sirupus ipecacuanhae; Sirupus plantaginis lanceolatae; Sirupus rubi idaei; Sirupus simplex; Sorbitol solution 70% (non-crystallizing); Sorbitolum liquidum non cristallisabile; Sugar syrup DAB

Special notes

When selecting the *Syrups* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Syrups*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Marshmallow syru...	12067103	30753	40	20120312

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Sorbitol solutio...	12124910	30870	40	20120525
Caelo	Ribwort plantain...	11380011	30897	60	20120112
Caelo	Raspberry syrup ...	13122003	31051	40	1402567
Caelo	Sorbitol solutio...	13060206	31124	40	20130306
Caelo	Marshmallow syru...	12303106	31147	40	1404662
Caelo	Cherry syrup (co...	13278301	31155	40	1402568
apomix	Sugar syrup DAB	73A101013	31175	40	1403164
Fagron	Ipecacuanha syru...	13G23-M02	31181	40	1404013
Caelo	Ribwort plantain...	13152307	31218	40	1404663
Caelo	Sorbitol solutio...	14034706	31353	40	8780
Caelo	Sugar syrup DAB	14134301	31355	40	20140514
Caelo	Cherry syrup (co...	14034503	31359	40	20140131
Caelo	Marshmallow syru...	13410907	31365	40	20140318
Caelo	Raspberry syrup ...	14114303	31467	40	20140411
Euro OTC	Raspberry syrup ...	1508015	31773	40	20150814*
Euro OTC	Raspberry syrup ...	1511024	31802	40	20150511*
Caelo	Sugar syrup DAB	153348	31877	40	20151014*
Euro OTC	Raspberry syrup ...	1605008	32284	40	20160513*
Euro OTC	Raspberry syrup ...	1701022	32984	40	20170126*
Caelo	Cherry syrup (co...	163453	33099	40	20161221*
Fagron	Ipecacuanha syru...	02B200616	33263	40	20160627
Caelo	Marshmallow syru...	171283	33331	40	20170511*
Fagron	Ipecacuanha syru...	01A070817	33642	40	20170810

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 980 spectra of 24 reference samples from the substance/substance group *Syrups*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 24 different batches.
- 34 564 spectra from a total of 816 batches from further 205 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 1650 spectra from 45 reference samples of the substance/substance group *Syrups*.
- These include spectra of independent samples from 36 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
zu erfassen	Raspberry syrup (conserved)	zu erfassen	30176	30
Caelo	Cherry syrup (conserved)	10194902	30235	40
apomix	Sugar syrup DAB	52A150711	30401	40
apomix	Sugar syrup DAB	12A010212	30443	40
Caelo	Cherry syrup (conserved)	11168902	30478	40
Fagron	Ipecacuanha syrup (NRF 19.1)	11J20-M02	30484	40
Caelo	Sorbitol solution 70% (non-c...)	14034706	31353 [†]	20
Caelo	Sugar syrup DAB	14134301	31355 [†]	20
Caelo	Cherry syrup (conserved)	14034503	31359 [†]	20
Caelo	Raspberry syrup (conserved)	14114303	31467 [†]	20
Euro OTC	Raspberry syrup (conserved)	1508015	31773 [†]	20
Euro OTC	Raspberry syrup (conserved)	1511024	31802 [†]	20
Caelo	Sugar syrup DAB	153348	31877 [†]	20
Euro OTC	Raspberry syrup (conserved)	1603007	32193	40
Caelo	Marshmallow syrup (conserved)	161233	32210	40
Caelo	Raspberry syrup (conserved)	161516	32213	40
Caelo	Marshmallow syrup (conserved)	160775	32215	40
Caelo	Sorbitol solution 70% (non-c...)	160788	32217	40
Euro OTC	Sorbitol solution 70% (non-c...)	1604020	32277	40
Caelo	Marshmallow syrup (conserved)	162136	32335	40
Euro OTC	Raspberry syrup (conserved)	1608003	32705	40
Euro OTC	Raspberry syrup (conserved)	1610026	32833	40
Caelo	Raspberry syrup (conserved)	163454	33092	40
Caelo	Sugar syrup DAB	170581	33126	40
Euro OTC	Sorbitol solution 70% (non-c...)	1702011	33228	40
Caelo	Raspberry syrup (conserved)	171284	33332	40
apomix	Ipecacuanha syrup (NRF 19.1)	03A171016	33420	40
apomix	Ipecacuanha syrup (NRF 19.1)	02A200616	33421	40
Euro OTC	Raspberry syrup (conserved)	1705037	33450	40
Caelo	Ribwort plantain syrup (cons...)	16335301	33457	40
Caelo	Ribwort plantain syrup (cons...)	16335301	33458	40
apomix	Sugar syrup DAB	30A130617	33505	40
Euro OTC	Raspberry syrup (conserved)	1709003	33660	40
Caelo	Cherry syrup (conserved)	16345303	33837	40
Caelo	Cherry syrup (conserved)	173511	33891	40
Caelo	Sorbitol solution 70% (non-c...)	173054	33895	40
Euro OTC	Raspberry syrup (conserved)	1712004	33930	40
Euro OTC	Sorbitol solution 70% (non-c...)	1802020	34006	40
Caelo	Ribwort plantain syrup (cons...)	18078901	34048	40
Euro OTC	Raspberry syrup (conserved)	1803009	34092	40
Euro OTC	Sorbitol solution 70% (non-c...)	1804019	34182	40
Euro OTC	Raspberry syrup (conserved)	1807026	34190	40
Caelo	Raspberry syrup (conserved)	181909	34222	40
Caelo	Sugar syrup DAB	181757	34246	40
Caelo	Sugar syrup DAB	182278	34247	40

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 92 224 spectra from a total of 2124 batches from further 430 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 918 spectra from 185 *Apo-Ident* customers from 459 batches from the substance/substance group *Syrups*.
- These include spectra of independent samples from 445 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Marshmallow syrup (conserved)	13152307	2
ACM	Marshmallow syrup (conserved)	1015	1
Caelo	Marshmallow syrup (conserved)	13410909	1
Caelo	Marshmallow syrup (conserved)	13212910	1
Caelo	Marshmallow syrup (conserved)	14102602	2
Caelo	Marshmallow syrup (conserved)	14364804	1
Caelo	Marshmallow syrup (conserved)	13212908	1
Caelo	Marshmallow syrup (conserved)	15418604	2
Caelo	Marshmallow syrup (conserved)	1548604	1
Caelo	Marshmallow syrup (conserved)	12303101	1
Caelo	Marshmallow syrup (conserved)	15418605	1
Caelo	Marshmallow syrup (conserved)	16077502	1
Caelo	Marshmallow syrup (conserved)	13212909	2
Caelo	Marshmallow syrup (conserved)	12303106	3
Sanacorp	Marshmallow syrup (conserved)	13212909	1
Caelo	Marshmallow syrup (conserved)	13410907	1
Caelo	Raspberry syrup (conserved)	5021510	1
Caelo	Raspberry syrup (conserved)	4101405	1
Euro OTC	Raspberry syrup (conserved)	1306006-01	3
Caelo	Raspberry syrup (conserved)	3121404	1
Euro OTC	Raspberry syrup (conserved)	11081402	1
Caelo	Raspberry syrup (conserved)	13051501	1
Euro OTC	Raspberry syrup (conserved)	1310016-01	10
Caelo	Raspberry syrup (conserved)	13338001	10
Caelo	Raspberry syrup (conserved)	13393601	3
Caelo	Raspberry syrup (conserved)	13393602	4
Euro OTC	Raspberry syrup (conserved)	140300201	1
noweda	Raspberry syrup (conserved)	13393605	1
Euro OTC	Raspberry syrup (conserved)	1406031-01	5
Phönix	Raspberry syrup (conserved)	14114303	2
Caelo	Raspberry syrup (conserved)	140903901	1
Caelo	Raspberry syrup (conserved)	14114303	1
Caelo	Raspberry syrup (conserved)	14114304	5
Eurho/Phönix	Raspberry syrup (conserved)	1409039-01	1
Caelo	Raspberry syrup (conserved)	14114305	3
phoenix 14.11.14	Raspberry syrup (conserved)	14114306	1
Caelo	Raspberry syrup (conserved)	14114306	2
Caelo	Raspberry syrup (conserved)	14243405	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Raspberry syrup (conserved)	14243403	1
Caelo	Raspberry syrup (conserved)	14243401	1
Caelo	Raspberry syrup (conserved)	14371104	1
Caelo	Raspberry syrup (conserved)	14371108	1
Caelo	Raspberry syrup (conserved)	14371106	1
Fagron	Raspberry syrup (conserved)	1501011-01	1
Noweda 04.08.15	Raspberry syrup (conserved)	1501011-01	1
Euro OTC	Raspberry syrup (conserved)	1504011-01	3
Euro OTC	Raspberry syrup (conserved)	1508015-01	3
Euro OTC	Raspberry syrup (conserved)	1511023-01	3
Euro OTC	Raspberry syrup (conserved)	1511024-01	4
Caelo	Raspberry syrup (conserved)	151102401	1
Caelo	Raspberry syrup (conserved)	15149201	2
Caelo	Raspberry syrup (conserved)	15149202	2
Caelo	Raspberry syrup (conserved)	15149202	3
Caelo	Raspberry syrup (conserved)	15149205	8
Caesar & Loretz	Raspberry syrup (conserved)	15149203	1
Caelo	Raspberry syrup (conserved)	15320801	1
Caelo	Raspberry syrup (conserved)	15320803	3
Caesar&Loretz	Raspberry syrup (conserved)	15320801	1
Caesar&Loretz	Raspberry syrup (conserved)	15320804	1
Caelo	Raspberry syrup (conserved)	15320804	1
Caelo	Raspberry syrup (conserved)	15320813	1
Caelo	Raspberry syrup (conserved)	16003601	4
Caesar&Loretz	Raspberry syrup (conserved)	16003602	1
Euro OTC	Raspberry syrup (conserved)	1603007-01	6
Euro OTC	Raspberry syrup (conserved)	1605008-01	8
EuRho	Raspberry syrup (conserved)	1603007-01	1
Euro OTC	Raspberry syrup (conserved)	160508-01	1
Euro OTC	Raspberry syrup (conserved)	160500801	1
Caelo	Raspberry syrup (conserved)	1608003-01	1
Euro OTC	Raspberry syrup (conserved)	1608003-01	1
Euro OTC	Raspberry syrup (conserved)	1610026-01	2
Caelo	Raspberry syrup (conserved)	161002601	1
Caelo	Raspberry syrup (conserved)	16151601	1
Caesar & Loretz/Noweda	Raspberry syrup (conserved)	16151603	1
Caelo	Raspberry syrup (conserved)	16345401	1
Caesar & Loretz	Raspberry syrup (conserved)	16151602	2
Caesar & Loretz	Raspberry syrup (conserved)	16345403	3
Euro OTC	Raspberry syrup (conserved)	1712004-01	1
Caelo	Raspberry syrup (conserved)	17128404	2
Euro OTC	Raspberry syrup (conserved)	1705037-01	1
Euro OTC	Raspberry syrup (conserved)	30039701/001	1
Caesar & Loretz/ Noweda	Raspberry syrup (conserved)	18096702	1
Caesar & Loretz/ Caesar...	Raspberry syrup (conserved)	17128409	1
Caesar & Loretz/ Caesar...	Raspberry syrup (conserved)	18080601	1
Caelo	Raspberry syrup (conserved)	13194301	7
noweda	Raspberry syrup (conserved)		1
Caelo	Raspberry syrup (conserved)	11021405	1
Caelo	Raspberry syrup (conserved)	1081201	1
Caelo	Raspberry syrup (conserved)	4071203	1
Caelo	Raspberry syrup (conserved)	1105022-01	1
Caelo	Raspberry syrup (conserved)	11016508	1
Eurho / Sanacorp	Raspberry syrup (conserved)	1110025-01	2
Caelo	Raspberry syrup (conserved)	3011301	1
Caelo	Raspberry syrup (conserved)	11169003	1
Caelo	Raspberry syrup (conserved)	8021301	1
Caelo	Raspberry syrup (conserved)	11291404	1
Caelo	Raspberry syrup (conserved)	12038501	3
Caelo	Raspberry syrup (conserved)	11363004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
phönix	Raspberry syrup (conserved)	12038501	1
Caelo	Raspberry syrup (conserved)	12038502	2
Caelo	Raspberry syrup (conserved)	12038505	1
Euro OTC	Raspberry syrup (conserved)	1205001-01	4
EuRho	Raspberry syrup (conserved)	1205001-01	3
Anzag 31.10.12	Raspberry syrup (conserved)	1205001-01	1
Euro OTC	Raspberry syrup (conserved)	1210004-01	14
EuRho Anzag	Raspberry syrup (conserved)	1210004-01	1
Noweda	Raspberry syrup (conserved)	1210004-01	1
Caelo	Raspberry syrup (conserved)	1210004-01	1
Caelo	Raspberry syrup (conserved)	12121203	1
Caelo	Raspberry syrup (conserved)	12161705	4
Caelo	Raspberry syrup (conserved)	12161703	1
Anzag 13.11.12	Raspberry syrup (conserved)	12283101	1
Caelo	Raspberry syrup (conserved)	12283101	3
Caelo	Raspberry syrup (conserved)	12283102	1
Sanacorp	Raspberry syrup (conserved)	12283103	3
Caelo	Raspberry syrup (conserved)	12283103	4
Gehe	Raspberry syrup (conserved)	12367602	1
Caelo	Raspberry syrup (conserved)	12367602	9
Caelo	Raspberry syrup (conserved)	12367601	1
Euro OTC	Raspberry syrup (conserved)	1302023-01	2
jenne eurho	Raspberry syrup (conserved)	1302026-01	1
EuRho/Sanacorp	Raspberry syrup (conserved)	1302026-01	1
Eurho / Sanacorp	Raspberry syrup (conserved)	1302026-01	2
Caelo	Raspberry syrup (conserved)	1302026-01	3
Euro OTC	Raspberry syrup (conserved)	1302026-01	1
Anzag	Raspberry syrup (conserved)	1306006-01	1
Caelo	Raspberry syrup (conserved)	13122001	10
Caesar & Loretz/ Noweda	Raspberry syrup (conserved)	13122003	1
Caelo	Raspberry syrup (conserved)	13122002	1
Caelo	Raspberry syrup (conserved)	13122003	2
Caelo	Raspberry syrup (conserved)	13194302	1
Phoenix	Raspberry syrup (conserved)	13338001	1
Phönix	Raspberry syrup (conserved)	141114303	1
Euro OTC	Raspberry syrup (conserved)	13802025-01	1
Caelo	Raspberry syrup (conserved)	15051312	1
Caelo	Raspberry syrup (conserved)	2x13122001	1
Caelo	Raspberry syrup (conserved)	3961/V07/14	1
Caelo	Raspberry syrup (conserved)	17101202	1
Caelo	Raspberry syrup (conserved)	24081201	1
Caelo	Raspberry syrup (conserved)	24091302	1
Caelo	Raspberry syrup (conserved)	15071302	1
Caelo	Raspberry syrup (conserved)	18111308	1
Apomix/Gehe	Ipecacuanha syrup (NRF 19.1)	3A170915	1
Caelo	Ipecacuanha syrup (NRF 19.1)	13G23-M02	1
Fagron	Ipecacuanha syrup (NRF 19.1)	13G23-M02	1
Fagron	Ipecacuanha syrup (NRF 19.1)	14D17-M01	1
Fagron	Ipecacuanha syrup (NRF 19.1)	15C11-M08	3
Fagron	Ipecacuanha syrup (NRF 19.1)	11B10-M01	1
Caelo	Cherry syrup (conserved)	13278301	2
Caelo	Cherry syrup (conserved)	14034502	2
Caelo	Cherry syrup (conserved)	14034503	1
Phönix	Cherry syrup (conserved)	14114304	1
Caelo	Cherry syrup (conserved)	14371001	1
Caelo	Cherry syrup (conserved)	16029301	2
Caelo	Cherry syrup (conserved)	15149402	1
Caelo	Cherry syrup (conserved)	16345301	2
Caelo	Cherry syrup (conserved)	16195502	1
Caelo	Cherry syrup (conserved)	12367501	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cherry syrup (conserved)	109212258	1
Anzag	Cherry syrup (conserved)	11291603	1
Caelo	Cherry syrup (conserved)	11168902	1
Caelo	Cherry syrup (conserved)	13166201	1
Caelo	Raspberry syrup (conserved)	14114302	1
Sanacorp	Raspberry syrup (conserved)	14243405	1
Caelo	Raspberry syrup (conserved)	14371105	1
Caelo	Raspberry syrup (conserved)	14371107	1
Euro OTC	Raspberry syrup (conserved)	121000401	1
Gehe	Raspberry syrup (conserved)	12161702	1
A.H.D.	Raspberry syrup (conserved)	12367602	2
AHD	Raspberry syrup (conserved)	13122001	1
Euro OTC	Raspberry syrup (conserved)	1302026-01-335031904	1
Sanacorp	Raspberry syrup (conserved)	13194301	1
Apomix/Gehe	Sugar syrup DAB	9A030214	2
apomix/Phönix	Sugar syrup DAB	2A130116	1
Caelo	Sugar syrup DAB	15427102	2
apomix/Phönix	Sugar syrup DAB	9A280116	1
Apomix/Anzag	Sugar syrup DAB	1314I-02263	1
Apomix/Gehe	Sugar syrup DAB	15A100215	1
Caelo	Sugar syrup DAB	16091501	2
Caelo	Sugar syrup DAB	16214002	2
Caelo	Sugar syrup DAB	16214101	2
Caelo	Sugar syrup DAB	17028002	4
Apomix	Sugar syrup DAB	21A050314	3
Gutenberg-Apotheke	Sugar syrup DAB	18A030314	1
Apomix/Gehe	Sugar syrup DAB	49A020715	2
Apomix/Gehe	Sugar syrup DAB	60A020914	1
Apomix/Gehe	Sugar syrup DAB	32A240415	1
Apomix/Gehe	Sugar syrup DAB	81A251114	1
Apomix/Gehe	Sugar syrup DAB	09A030214	1
Phoenix/Apomix	Sugar syrup DAB	29A220313	1
apomix/Phoenix	Sugar syrup DAB	20A280213	1
Apommix-Sanacorp	Sugar syrup DAB	51A090713-3399310141	1
Anzag	Sugar syrup DAB	44A260612	1
apomix/Phoenix	Sugar syrup DAB	73A101013	1
Apomix/Sanacorp	Sugar syrup DAB	66A110913-3442311102	1
Caelo	Sorbitol solution 70% (non-c...	3031404	1
Caelo	Sorbitol solution 70% (non-c...	12/18-N09	1
Caelo	Sorbitol solution 70% (non-c...	1207012-02	2
Euro OTC	Sorbitol solution 70% (non-c...	1207012-02	3
Caelo	Sorbitol solution 70% (non-c...	13060204	5
Euro OTC	Sorbitol solution 70% (non-c...	1309026-01	5
Ichthyol/Gehe	Sorbitol solution 70% (non-c...	1309026-01	1
Caelo	Sorbitol solution 70% (non-c...	1309026-01	2
Euro OTC	Sorbitol solution 70% (non-c...	1309026-02	5
EuRho	Sorbitol solution 70% (non-c...	1309026-02	1
Phoenix	Sorbitol solution 70% (non-c...	1309026-02	1
Caelo	Sorbitol solution 70% (non-c...	13233801	2
Caelo	Sorbitol solution 70% (non-c...	13233804	1
Caelo	Sorbitol solution 70% (non-c...	13233803	3
Caelo	Sorbitol solution 70% (non-c...	13233807	3
Caelo	Sorbitol solution 70% (non-c...	13233808	3
Caelo	Sorbitol solution 70% (non-c...	13233810	6
Gehe	Sorbitol solution 70% (non-c...	13233810	1
Fagron	Sorbitol solution 70% (non-c...	13E24-B01	2
Fagron	Sorbitol solution 70% (non-c...	13J04-B01-289676	1
Caelo	Sorbitol solution 70% (non-c...	14034702	1
Caelo	Sorbitol solution 70% (non-c...	14034706	1
Euro OTC	Sorbitol solution 70% (non-c...	1412031-01	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Sorbitol solution 70% (non-c...	141203101	1
Caelo	Sorbitol solution 70% (non-c...	1412031-01	1
Euro OTC	Sorbitol solution 70% (non-c...	1412031-02	4
Caelo	Sorbitol solution 70% (non-c...	14169401	1
Caelo	Sorbitol solution 70% (non-c...	14169403	4
Caelo	Sorbitol solution 70% (non-c...	14169409	2
Caelo	Sorbitol solution 70% (non-c...	14298101	7
Caelo	Sorbitol solution 70% (non-c...	14298102	3
Caelo	Sorbitol solution 70% (non-c...	15159301	5
Caelo	Sorbitol solution 70% (non-c...	15159302	6
Caelo	Sorbitol solution 70% (non-c...	15175301	1
Caelo	Sorbitol solution 70% (non-c...	15175306	3
Caelo	Sorbitol solution 70% (non-c...	15175307	1
Caelo	Sorbitol solution 70% (non-c...	15175309	4
Caelo	Sorbitol solution 70% (non-c...	1519302	1
Caelo	Sorbitol solution 70% (non-c...	15175310	3
Caelo	Sorbitol solution 70% (non-c...	15437900	1
Caelo	Sorbitol solution 70% (non-c...	15437906	4
Caelo	Sorbitol solution 70% (non-c...	15437907	2
Euro OTC	Sorbitol solution 70% (non-c...	1604020-01	7
Euro OTC	Sorbitol solution 70% (non-c...	1604020-02	1
Euro OTC	Sorbitol solution 70% (non-c...	1604020-03	1
Caelo	Sorbitol solution 70% (non-c...	16078806	3
Caelo	Sorbitol solution 70% (non-c...	16078807	1
Caelo	Sorbitol solution 70% (non-c...	16078810	2
Fagron	Sorbitol solution 70% (non-c...	16078810	1
Caelo	Sorbitol solution 70% (non-c...	16078811	1
Fagron	Sorbitol solution 70% (non-c...	29081505	1
Caelo	Sorbitol solution 70% (non-c...	17027010	2
Caelo	Sorbitol solution 70% (non-c...	17027005	1
Fagron	Sorbitol solution 70% (non-c...	21081406	1
Caelo	Sorbitol solution 70% (non-c...	17305401	1
Fagron	Sorbitol solution 70% (non-c...	2013M-02299	1
Caelo	Sorbitol solution 70% (non-c...	2071305	1
Caelo	Sorbitol solution 70% (non-c...	8051306	1
Caelo	Sorbitol solution 70% (non-c...	2071306	1
Gehe	Sorbitol solution 70% (non-c...	10323008	1
Euro OTC	Sorbitol solution 70% (non-c...	1102001-02	1
Caelo	Sorbitol solution 70% (non-c...	1102001-02	1
Caelo	Sorbitol solution 70% (non-c...	11255613	1
Fagron	Sorbitol solution 70% (non-c...	11G18-N14	1
Euro OTC	Sorbitol solution 70% (non-c...	1207012-01	4
Euro OTC	Sorbitol solution 70% (non-c...	120701201	1
Anzag	Sorbitol solution 70% (non-c...	120701201	1
Caelo	Sorbitol solution 70% (non-c...	1207012-01	2
Gehe	Sorbitol solution 70% (non-c...	1207012-01	1
EuRho	Sorbitol solution 70% (non-c...	1207012-02	1
Gehe	Sorbitol solution 70% (non-c...	1207012-02	1
Gehe	Sorbitol solution 70% (non-c...	120701202	1
Caelo	Sorbitol solution 70% (non-c...	12124906	2
Caelo	Sorbitol solution 70% (non-c...	12124905	1
Caelo	Sorbitol solution 70% (non-c...	12124903	2
Caelo	Sorbitol solution 70% (non-c...	12124908	2
Gehe	Sorbitol solution 70% (non-c...	12124910	1
Caelo	Sorbitol solution 70% (non-c...	12124912	4
Caelo	Sorbitol solution 70% (non-c...	12124914	5
Caelo	Sorbitol solution 70% (non-c...	12124915	1
Fagron	Sorbitol solution 70% (non-c...	12E08-N09	1
Gehe	Sorbitol solution 70% (non-c...	12E08-N09	1
Gehe	Sorbitol solution 70% (non-c...	12K07-N04	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Sorbitol solution 70% (non-c...	12118-N09	1
Caelo	Sorbitol solution 70% (non-c...	13060203	6
Caesar u Loretz	Sorbitol solution 70% (non-c...	13060206	1
Caelo	Sorbitol solution 70% (non-c...	13060206	2
Caelo	Sorbitol solution 70% (non-c...	13060211	1
Fagron	Sorbitol solution 70% (non-c...	139026	1
Caelo	Sorbitol solution 70% (non-c...	22011306	1
Caelo	Sorbitol solution 70% (non-c...	28081301	1
Caesar & Loretz GmbH	Sorbitol solution 70% (non-c...	14298101	1
Caelo	Sorbitol solution 70% (non-c...	26071209	1
EuRho/Phönix	Sorbitol solution 70% (non-c...	1412031-02	1
Caelo	Ribwort plantain syrup (cons...	13152304	1
Caelo	Ribwort plantain syrup (cons...	13152308	3
Caelo	Ribwort plantain syrup (cons...	14125702	2
Caelo	Ribwort plantain syrup (cons...	15076003	1
Sanacorp	Ribwort plantain syrup (cons...	13152303	1
Caelo	Ribwort plantain syrup (cons...	13152305	1
Caelo	Ribwort plantain syrup (cons...	12347501	1
Caelo	Ribwort plantain syrup (cons...	12347506	2
Caelo	Ribwort plantain syrup (cons...	13152302	1
apomix	Sugar syrup DAB	2A130116	1
Phönix Gotha	Sugar syrup DAB	1A130114	1
Apomix/Alliance Healthc...	Sugar syrup DAB	1A130116	1
apomix/alliance healthc...	Sugar syrup DAB	3a150116	1
Apomix/Kehr	Sugar syrup DAB	3A150116	1
apomix/alliance healthc...	Sugar syrup DAB	3A150116	1
Apomix/Kehr	Sugar syrup DAB	3A180118	5
AMH Niemann/Phoenix	Sugar syrup DAB	3A270114	1
Apomix/ Noweda	Sugar syrup DAB	5A220115	1
Apomix / Kehr	Sugar syrup DAB	5A220115	1
Apomix/Sanacorp	Sugar syrup DAB	2A130114	1
Apomix/Alliance Healthc...	Sugar syrup DAB	7A230115	1
Apomix / Kehr	Sugar syrup DAB	7A230117	1
Apomix/Sanacorp	Sugar syrup DAB	7A230115	1
AMH Niemann GmbH & Co.K...	Sugar syrup DAB	7A250116	1
Apomix PKH/V.d.L.	Sugar syrup DAB	7A310114	1
Apomix/Kehr	Sugar syrup DAB	8A070218	2
APOMIX/ Alliance Health...	Sugar syrup DAB	7A310114	1
Caelo	Sugar syrup DAB	9A280116	3
Noweda/Apomix	Sugar syrup DAB	9A280116	1
Apomix/Noweda	Sugar syrup DAB	9A250117	1
Apomix/Anzag	Sugar syrup DAB	1116E-02263	1
PKH-Halle/Sanacorp	Sugar syrup DAB	10A030214	1
Apomix / Kehr	Sugar syrup DAB	11A120216	2
apomix/Kehr	Sugar syrup DAB	11A280115	2
Alliance health	Sugar syrup DAB	11A280115	1
apomix	Sugar syrup DAB	13101405	1
Euro OTC	Sugar syrup DAB	13A150216	1
Klenk	Sugar syrup DAB	13A150216	2
AMH Niemann/Phoenix	Sugar syrup DAB	12A050214	1
apomix PKH GmbH Halle	Sugar syrup DAB	13A290115	1
Caelo	Sugar syrup DAB	14299801	4
Caelo	Sugar syrup DAB	14189102	2
Caelo	Sugar syrup DAB	13A1502216	1
Caelo	Sugar syrup DAB	14189101	1
Caelo	Sugar syrup DAB	14246001	1
Sanacorp	Sugar syrup DAB	13A290115	1
Caelo	Sugar syrup DAB	14337004	1
Caelo	Sugar syrup DAB	15117002	1
Apomix/Anzag	Sugar syrup DAB	1515Q-02263	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sugar syrup DAB	15334802	1
Apomix/Alliance Healthc...	Sugar syrup DAB	14A150216	1
Caelo	Sugar syrup DAB	15315201	1
Caelo	Sugar syrup DAB	16187902	2
Caelo	Sugar syrup DAB	15150602	1
JENNE	Sugar syrup DAB	16214101	1
Caelo	Sugar syrup DAB	17129001	2
Caelo	Sugar syrup DAB	17304302	2
Apomix/GEHE	Sugar syrup DAB	17A030314	1
Caelo	Sugar syrup DAB	17200501	1
Caelo	Sugar syrup DAB	17A240417	1
apomix	Sugar syrup DAB	18081517	1
Caelo	Sugar syrup DAB	18032801	1
Apomix/Gehe	Sugar syrup DAB	18A030314	1
Phoenix/apomix	Sugar syrup DAB	19A040314	1
Apomix / Kehr	Sugar syrup DAB	18A030314	1
Apomix/Kehr	Sugar syrup DAB	19A150517	1
PKH Halle	Sugar syrup DAB	19A040314	1
apomix	Sugar syrup DAB	19A040314	1
Apomix/Anzag	Sugar syrup DAB	2112I-02263	1
Caelo	Sugar syrup DAB	21A170316	1
Apomix/Alliance Healthc...	Sugar syrup DAB	20A270215	2
apomix	Sugar syrup DAB	21A130315	1
AMH Niemann GmbH & Co.KG	Sugar syrup DAB	22A130315	1
Caelo	Sugar syrup DAB	21A050314	1
Apomix	Sugar syrup DAB	22A050314	1
Apomix/Alliance Healthc...	Sugar syrup DAB	22A170316	1
Caelo	Sugar syrup DAB	22A130315	2
Apomix/Kehr	Sugar syrup DAB	23A290517	2
AHD / apomix	Sugar syrup DAB	24A180316	1
PKH/Phoenix	Sugar syrup DAB	25A310517	1
Apomix/Alliance Healthc...	Sugar syrup DAB	25A170314	1
Apomix/Kehr	Sugar syrup DAB	27A090418	1
apomix/Phoenix	Sugar syrup DAB	25A080415	1
Caelo	Sugar syrup DAB	23A180316	1
Apomix / Kehr	Sugar syrup DAB	27A110416	3
rPhR Neidhar	Sugar syrup DAB	30A190416	1
Apomix/Kehr	Sugar syrup DAB	29A130617	2
Caelo	Sugar syrup DAB	30A190416	1
Apomix/Kehr	Sugar syrup DAB	30A200415	3
AMH Niemann GmbH	Sugar syrup DAB	30A200415	1
Sanacorp/ apomix	Sugar syrup DAB	28A110416	1
Apomix / Noweda	Sugar syrup DAB	30A190416	1
AMH Niemann GmbH & Co.K...	Sugar syrup DAB	33A170516	1
hochhaus-Apotheke	Sugar syrup DAB	35A300516	2
Caelo	Sugar syrup DAB	33A190617	1
Apomix/ Phönix	Sugar syrup DAB	32A240415	1
Apomix/Sanacorp	Sugar syrup DAB	35A070414	1
Kehr/Apomix	Sugar syrup DAB	35A310518	2
Alliance H.	Sugar syrup DAB	37/A030616	1
Apomix/Anzag	Sugar syrup DAB	37A03616	1
apomix	Sugar syrup DAB	37A220515	1
Apomix	Sugar syrup DAB	39A050614	1
apomix/Noweda	Sugar syrup DAB	39A050614	2
Caelo	Sugar syrup DAB	38A190514	1
Phönix Neidhardt	Sugar syrup DAB	40A050614	1
Apomix/ Phönix	Sugar syrup DAB	40A290515	2
apomix	Sugar syrup DAB	41A040716	1
Apomix/Kehr	Sugar syrup DAB	42A120615	1
Apomix/Kehr	Sugar syrup DAB	44A070817	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix	Sugar syrup DAB	41A120615	1
Caelo	Sugar syrup DAB	42A120615	1
Apomix / Noweda	Sugar syrup DAB	43A110716	2
Apomix/Alliance Healthc...	Sugar syrup DAB	44A110716	1
Apomix/Kehr	Sugar syrup DAB	45A210817	1
Apomix / Kehr	Sugar syrup DAB	48A210714	1
apomix	Sugar syrup DAB	49A220714	1
Apomix	Sugar syrup DAB	48A210714	1
Noweda / Apomix	Sugar syrup DAB	47A120613	1
apomix/Phoenix	Sugar syrup DAB	49A220714	2
PKH GmbH / Phoenix	Sugar syrup DAB	49A020715	1
Apomix / Kehr	Sugar syrup DAB	49A290716	3
AMH Niemann GmbH / Sana...	Sugar syrup DAB	49A290716	1
Apomix/Anzag	Sugar syrup DAB	50A290716	1
Apomix / Kehr	Sugar syrup DAB	50A020715	1
AMH Niemann GmbH	Sugar syrup DAB	50A020715	1
Apomix/Phönix	Sugar syrup DAB	51A070814	2
Apomix/Kehr	Sugar syrup DAB	51A010816	1
Apomix/Sancorp	Sugar syrup DAB	55A130814	1
Caelo	Sugar syrup DAB	51A070814	1
Apomix/Alliance Healthc...	Sugar syrup DAB	55A130814	1
apomix	Sugar syrup DAB	57A140814	1
Apomix	Sugar syrup DAB	55A130814	1
Apomix/Noweda	Sugar syrup DAB	57A180816	1
PKH Halle / Phoenix	Sugar syrup DAB	57A180816	1
apomix	Sugar syrup DAB	55A130814	1
Apomix	Sugar syrup DAB	55a220816	1
Apomix/Kehr	Sugar syrup DAB	57A260815	2
Apomix/Alliance Healthc...	Sugar syrup DAB	58A260815	1
Phoenix	Sugar syrup DAB	59A01914	1
apomix	Sugar syrup DAB	59A010914	1
Apomix/Phönix	Sugar syrup DAB	59A270815	1
Apomix/Kehr	Sugar syrup DAB	63A061117	3
apomix	Sugar syrup DAB	59A070814	1
Apomix / Kehr	Sugar syrup DAB	63A260916	4
Hochhaus-Apotheke	Sugar syrup DAB	64A260916	3
Apomix/Kehr	Sugar syrup DAB	66A141117	2
Apomix/Anzag	Sugar syrup DAB	67A131016	1
Apomix/Alliance Healthc...	Sugar syrup DAB	68A131016	1
Apomix / Kehr	Sugar syrup DAB	66A160914	2
hochhaus-Apotheke	Sugar syrup DAB	67A-131016	1
Apomix / Noweda	Sugar syrup DAB	67A131016	1
hochhaus-Apotheke	Sugar syrup DAB	67A131016	1
Apomix/ Phönix	Sugar syrup DAB	68A250915	1
Apomix / Kehr	Sugar syrup DAB	70328/4Hk	1
Apomix/Sanacorp	Sugar syrup DAB	72A041116	1
Apomix/Noweda	Sugar syrup DAB	71A291117	1
AMH Niemann GmbH	Sugar syrup DAB	70A051015	1
Apomix/Alliance Healthc...	Sugar syrup DAB	68A230914	1
Apomix/Alliance Healthc...	Sugar syrup DAB	71A091015	1
Apomix / Kehr	Sugar syrup DAB	72A041116	1
apomix/ Noweda	Sugar syrup DAB	72A091015	1
Apomix/Kehr	Sugar syrup DAB	72A091015	1
Apomix/Phönix	Sugar syrup DAB	72A091015	1
Phönix/Apomix	Sugar syrup DAB	73A101013	2
Apomix/Anzag	Sugar syrup DAB	76A231015	2
Kehr/Apomix	Sugar syrup DAB	74A241014	2
Apomix/Phönix	Sugar syrup DAB	74A151015	1
apomix Fiebig	Sugar syrup DAB	77A261015	1
Caelo	Sugar syrup DAB	77A061114	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix Gotha	Sugar syrup DAB	77A171013	1
Caelo	Sugar syrup DAB	79A281013	2
PKH / Phoenix	Sugar syrup DAB	80A061115	2
Caelo	Sugar syrup DAB	81A161115	4
Euro OTC	Sugar syrup DAB	80A061115	1
Apomix/Kehr	Sugar syrup DAB	85A261115	1
Caelo	Sugar syrup DAB	81A-251114	1
AMH Niemann GmbH	Sugar syrup DAB	85A261115	1
Caesar & Loretz GmbH	Sugar syrup DAB	86A051214	1
Apomix / Kehr	Sugar syrup DAB	86A051214	1
Apomix / Gehe	Sugar syrup DAB	87A071215	1
Kehr/Apomix	Sugar syrup DAB	88A121214	2
Apomix / Kehr	Sugar syrup DAB	88A12124	1
Noweda/Apomix	Sugar syrup DAB	87A051214	1
Apomix/Gehe	Sugar syrup DAB	89A111215	1
apomix VDL	Sugar syrup DAB	93A111213	1
apomix/ Phönix	Sugar syrup DAB	93A231214	1
Apomix / Noweda	Sugar syrup DAB	93A231214	1
apomix/Sanacorp	Sugar syrup DAB	96A161213	1
Phönix/apomix	Sugar syrup DAB	97A191213	2
Apomix	Sugar syrup DAB	1A140113	1
PKH / Phoenix	Sugar syrup DAB	93A231214	1
AMH Niemann	Sugar syrup DAB	1A130114	1
Apomix	Sugar syrup DAB	99A201213	1
Apomix PKH/V.d.L.	Sugar syrup DAB	9A290113	1
Apomix	Sugar syrup DAB	12A010212	1
Noweda	Sugar syrup DAB	11A060213	1
Apomix	Sugar syrup DAB	5A210113	1
apomix	Sugar syrup DAB	1A140115	1
Apomix	Sugar syrup DAB	10A3000113	1
Caelo	Sugar syrup DAB	13A120213	1
apomix/Phoenix	Sugar syrup DAB	15A210213	4
apomix/ Phoenix	Sugar syrup DAB	11A060213	1
apomix	Sugar syrup DAB	11121305	1
Fagron	Sugar syrup DAB	13089R	1
apomix	Sugar syrup DAB	16A210213	1
apomix/Phönix	Sugar syrup DAB	214050314	1
apomix	Sugar syrup DAB	20A270215	1
Phönix/Apomix	Sugar syrup DAB	21A010313	2
apomix	Sugar syrup DAB	20021307	1
apomix/ Phoenix	Sugar syrup DAB	21A010313	1
Apomix	Sugar syrup DAB	17A270213	1
Apomix/Sancorp	Sugar syrup DAB	21A010313	1
Apomix	Sugar syrup DAB	25A110313	1
Caelo	Sugar syrup DAB	21A010313	1
Apomix/Phönix	Sugar syrup DAB	26A190313	1
apomix	Sugar syrup DAB	23081202	1
Apomix	Sugar syrup DAB	25A200412	1
Sanacorp/ PKH Halle	Sugar syrup DAB	29A110512	1
apomix	Sugar syrup DAB	26A080415	1
apomix	Sugar syrup DAB	27A270712	1
apomix	Sugar syrup DAB	27A200313	1
apomix/Phoenix	Sugar syrup DAB	29A220313	3
Apomix	Sugar syrup DAB	32A120413	1
Apomix/Phönix	Sugar syrup DAB	33A120413	3
Apomix/ Phönix	Sugar syrup DAB	34A240413	1
apomix	Sugar syrup DAB	35A250512	1
Apomix PKH/V.d.L.	Sugar syrup DAB	37A240413	2
apomix	Sugar syrup DAB	36A140414	1
Fagron	Sugar syrup DAB	38A250413	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix	Sugar syrup DAB	43A040613	1
PKH	Sugar syrup DAB	41A190612	1
Apomix / Noweda	Sugar syrup DAB	41A230513	1
Apomix, Noweda	Sugar syrup DAB	43A040613	2
Nachprüfung	Sugar syrup DAB	42A220612	1
apomix	Sugar syrup DAB	43A260612	1
Caelo	Sugar syrup DAB	44A050613	1
apomix	Sugar syrup DAB	45A050613	1
apomix	Sugar syrup DAB	46A030712	1
Apomix	Sugar syrup DAB	47A050712	3
Apomix/Sanacorp	Sugar syrup DAB	45A050613	1
Apomix Anzag	Sugar syrup DAB	49A100712	1
Phönix Gotha	Sugar syrup DAB	51A090713	3
Gehe	Sugar syrup DAB	51A090713	1
Apomix / Noweda	Sugar syrup DAB	51A090713	1
Phönix Gotha	Sugar syrup DAB		1
Niemann GmbH & Co.KG	Sugar syrup DAB	53A080814	1
Apomix/ Gehe	Sugar syrup DAB	55A180712	2
apomix/Phoenix	Sugar syrup DAB	56A160813	3
noweda/Apomix	Sugar syrup DAB	56A160813	1
Sanacorp	Sugar syrup DAB	58A270712	1
Anzag	Sugar syrup DAB	58A210813	1
Phönix/Apomix	Sugar syrup DAB	58A270712	1
Gehe/ Apomix	Sugar syrup DAB	59A030812	1
Phönix/Lamotte	Sugar syrup DAB	60A230813	2
Phönix Gotha	Sugar syrup DAB	59A220813	1
apomix	Sugar syrup DAB	58A270712	1
Apomix PKH/V.d.L.	Sugar syrup DAB	64A050913	1
Caelo	Sugar syrup DAB	58A270712	1
apomix	Sugar syrup DAB	66A210912	1
Phönix/apomix	Sugar syrup DAB	67A110913	2
Phönix/ apomix	Sugar syrup DAB	66A110913	1
apomix	Sugar syrup DAB	67A240912	1
apomix	Sugar syrup DAB	71A091015	1
PKH-Halle/Sanacorp	Sugar syrup DAB	72A041013	2
Anzag 25.01.2013	Sugar syrup DAB	71A091012	1
AMH Niemann GmbH	Sugar syrup DAB	74A241014	1
Caelo	Sugar syrup DAB	76A171013	1
apomix	Sugar syrup DAB	76A061114	1
Apomix/Noweda	Sugar syrup DAB	78A011112	1
apomix/ Phoenix	Sugar syrup DAB	76A301012	1
Sanacorp	Sugar syrup DAB	82A081112	1
Phönix/apomix	Sugar syrup DAB	82A081112	3
Apomix	Sugar syrup DAB	81A-081112	1
apomix	Sugar syrup DAB	83A191115	1
apomix	Sugar syrup DAB	83A131112	1
Anzag 17.01.2013	Sugar syrup DAB	88A051212	1
Apomix /Alliance	Sugar syrup DAB	85A191113	1
Apomix	Sugar syrup DAB	93A111213	1
apomix/Phoenix	Sugar syrup DAB	89A271113	1
Apomix/Gehe	Sugar syrup DAB	93A271212	1
Apomix	Sugar syrup DAB	94A281212	1
apomix	Sugar syrup DAB	95A111213	1

- 90 434 spectra from 501 *Apo-Ident* customers from a total of 26 564 batches from a further 327 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Syrups* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Syrups* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	980	0	34 564
Type B	0	1640	10	92 224
Type C	0	907	11	90 434

The substance/substance group *Syrups* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.3878 %)
Type B	100.0000 % (> 99.9761 %)	99.3939 % (> 99.2121 %)
Type C	100.0000 % (> 99.6586 %)	98.8017 % (> 98.4749 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30753	31147	3.53	32.27
30870	31353	7.97	24.32
30897	31877	3.44	32.98
31051	31051	0.00	24.43
31124	31353	2.00	25.20
31147	31147	0.00	30.92
31155	31155	0.00	24.14
31175	31175	0.00	30.01
31181	31181	0.00	28.98
31218	31218	0.00	34.22
31353	31353	0.00	25.20
31355	31175	1.89	31.00
31359	31051	1.46	25.18

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31365	31877	3.23	32.17
31467	31051	1.83	24.63
31773	31773	0.00	27.84
31802	31802	0.00	25.78
31877	31877	0.00	31.73
32284	32284	0.00	25.42
32984	32984	0.00	28.38
33099	33099	0.00	22.73
33263	31181	5.72	30.27
33331	33331	0.00	26.93
33642	31181	2.48	27.32

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Tea tree oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31023-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Tea tree oil; Melaleuca alternifolia oleum; Oleum melaleuca alternifolia

Special notes

When selecting the *Tea tree oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	European Pharmacopoeia 8th Edition, Basic Version 2014 [3]
Komm2.2.40	Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident [4]
AA004	Erstellung und Validierung eines IdentModul-Updates
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Tea tree oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Tea tree oil	12201303	31023	40	1404597
Caelo	Tea tree oil	13288914	31401	40	20130828
Caelo	Tea tree oil	14203203	31960	40	20140703

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Tea tree oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 179 spectra from 5 reference samples of the substance/substance group *Tea tree oil*.
- These include spectra of independent samples from 4 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Tea tree oil	00244265	30306	39
Caelo	Tea tree oil	14203203	31960 [†]	20
Caelo	Tea tree oil	15337105	33478	40
Caelo	Tea tree oil	173227	34058	40
Caelo	Tea tree oil	170677	34059	40

- 93 695 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 86 spectra from 31 *Apo-Ident* customers from 45 batches from the substance/substance group *Tea tree oil*.
- These include spectra of independent samples from 40 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Tea tree oil	11333516	1
Caelo	Tea tree oil	14203202	1
primavera	Tea tree oil	5604136	1
Taoasis	Tea tree oil	41015006139	1
Caelo	Tea tree oil	11358910	1
Caelo	Tea tree oil	12201303	1
Caelo	Tea tree oil	12201303	2
Caelo	Tea tree oil	12201304	3
Caelo	Tea tree oil	12201311	4
Caelo	Tea tree oil	12201313	1
Apoth.Bauer&Co,Münch	Tea tree oil		2
Caelo	Tea tree oil	13099506	6
Caelo	Tea tree oil	13099510	4
Caelo	Tea tree oil	13288902	1
Caelo	Tea tree oil	14071604	2
Caelo	Tea tree oil	14111201	1
Caelo	Tea tree oil	14203203	4
Caelo	Tea tree oil	14302301	8
Caelo	Tea tree oil	14203201	1
Caelo	Tea tree oil	15051806	2
Caelo	Tea tree oil	14302302	1
Caelo	Tea tree oil	15051807	2
Caelo	Tea tree oil	15337106	3
Minck	Tea tree oil	1511E-03093	1
Ebert 01.09.2016	Tea tree oil	15337105	1
Caelo	Tea tree oil	15337111	2
Caelo	Tea tree oil	15337112	4
Caelo	Tea tree oil	15337120	1
Caelo	Tea tree oil	170667706	2
Bombastus	Tea tree oil	304033	1
Taoasis	Tea tree oil	516495-108562	1
Caelo	Tea tree oil	16007506	1
Taoasis	Tea tree oil	33152-112369	1
Taoasis	Tea tree oil	535627-110746	2
Taoasis	Tea tree oil	556031-114257	2
Fagron	Tea tree oil		2
Apotheker bauer	Tea tree oil		1
apotheker bauer& cie oHG	Tea tree oil	6548331306	1
Apoth.Bauer&Co,Münch	Tea tree oil	692624	1
Bergland-Pharma/Anzag	Tea tree oil	C002685	1
Bergland-Pharma/Anzag	Tea tree oil	C101782	1
Sensient	Tea tree oil	C9663	1
Taoasis GmbH	Tea tree oil	98546	1
Gehe	Tea tree oil	K203705	3
Bergland /Noweda	Tea tree oil	K302332	1
Bergland / Noweda	Tea tree oil	P102633	1

- 91 266 spectra from 501 *Apo-Ident* customers from a total of 26 965 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Tea tree oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Tea tree oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	170	9	93 695
Type C	0	77	9	91 266

The substance/substance group *Tea tree oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	94.9721 % (> 93.2961 %)
Type C	100.0000 % (> 99.6587 %)	89.5349 % (> 86.0465 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31023	31023	0.00	50.69
31401	31023	25.58	48.52
31960	31023	7.65	52.04

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Thyme oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30206-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Thyme oil; Oleum thymi

Special notes

When selecting the *Thyme oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Thyme oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Thyme oil	12020518	31018	40	1407220
Caelo	Thyme oil	14225405	31566	40	20140822
Caelo	Thyme oil	15230104	32265	40	20150702
Caelo	Thyme oil	163279	33083	40	20161222*
Bombastus	Thyme oil	17001282	33878	40	20171108

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 200 spectra of 5 reference samples from the substance/substance group *Thyme oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 5 different batches.
- 35 344 spectra from a total of 835 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 140 spectra from 4 reference samples of the substance/substance group *Thyme oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Thyme oil	10078906	30206	40
Caelo	Thyme oil	14225405	31566 [†]	20
Caelo	Thyme oil	16327902	33477	40
Bombastus	Thyme oil	305865	33877	40

- 93 734 spectra from a total of 2163 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 155 spectra from 76 *Apo-Ident* customers from 65 batches from the substance/substance group *Thyme oil*.
- These include spectra of independent samples from 61 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Thyme oil	13044710	1
Caelo	Thyme oil	13044716	3
VDL	Thyme oil	14225413	1
Caelo	Thyme oil	14225414	1
Caelo	Thyme oil	14225403	1
Caelo	Thyme oil	154022002	1
Caelo	Thyme oil	15402204	1
Caelo	Thyme oil	15402202	16
Caelo	Thyme oil	14225418	1
Caelo	Thyme oil	12020515	1
Caelo	Thyme oil	12020518	7
Caelo	Thyme oil	3313U-02023	1
Caelo	Thyme oil	13044705	1
Caelo	Thyme oil	12020514	5
Bombastus	Thyme oil	283494	1
Bombastus	Thyme oil	268010	2
Caelo	Thyme oil	13044707	4
Caelo	Thyme oil	6081401	1
PVL	Thyme oil	271112	1
Taoasis	Thyme oil	109378	1
Caelo	Thyme oil	10118603	1
Phönix 03.12.12	Thyme oil	11275811	1
Caelo	Thyme oil	11275811	1
Caelo	Thyme oil	11275806	1
Caelo	Thyme oil	12020503	4
Fiebi 07.11.12 ek:8,48	Thyme oil	12020503	1
Caelo	Thyme oil	11275801	1
Fagron	Thyme oil	10L20-N09	1
Anzag 27.12.12 ek:8,71	Thyme oil	12020510	2
Caelo	Thyme oil	12020510	7
Caelo	Thyme oil	12020511	1
Caelo	Thyme oil	1202510	1
Caelo	Thyme oil	12165401	1
Großhandel	Thyme oil	12165401	1
Caelo	Thyme oil	12201104	5
Caelo	Thyme oil	13021303	1
Caelo	Thyme oil	13044702	3
Caelo	Thyme oil	1221104	1
Caelo	Thyme oil	13044708	6
Bombastus	Thyme oil	13044708	1
	Thyme oil	13044708	1
02.07.2014 Phoenix	Thyme oil	13044715	3
Caelo	Thyme oil	13044709	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Thyme oil	13044715	2
Caelo	Thyme oil	14071202	1
Caelo	Thyme oil	142225412	1
Caelo	Thyme oil	14225401	3
Caelo	Thyme oil	14225405	2
Caelo	Thyme oil	14225408	1
Caelo	Thyme oil	14225409	2
Caelo	Thyme oil	14225406	5
Caelo	Thyme oil	14225413	5
Caelo	Thyme oil	14225410	2
Caelo	Thyme oil	14225412	2
Caelo	Thyme oil	14225417	3
Caelo	Thyme oil	15230104	6
Caelo	Thyme oil	1540220	1
Caesar & Loretz/Noweda	Thyme oil	15402202	1
Caelo	Thyme oil	15402203	1
Caelo	Thyme oil	15402207	3
Caelo	Thyme oil	15402208	2
Caelo	Thyme oil	16327906	1
Caelo	Thyme oil	15402219	1
Caelo	Thyme oil	16327909	1
Caelo	Thyme oil	16327910	2
Caelo	Thyme oil	2615A-02023	1
Post-Apotheke	Thyme oil	268011	2
Taoasis	Thyme oil	28739-104426	1
Taoasis	Thyme oil	28739-106640	1
Fagron	Thyme oil	606549	1
Taoasis	Thyme oil	31911-109378	1
Caelo	Thyme oil	3211I-02023	1
Taoasis	Thyme oil	31911-110649	1
Euro OTC	Thyme oil		1

- 91 197 spectra from 501 *Apo-Ident* customers from a total of 26 945 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Thyme oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Thyme oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	200	0	35 344
Type B	0	100	40	93 734
Type C	0	131	24	91 197

The substance/substance group *Thyme oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.0000 %)
Type B	100.0000 % (> 99.9761 %)	71.4286 % (> 69.2857 %)
Type C	100.0000 % (> 99.6586 %)	84.5161 % (> 82.5806 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31018	31018	0.00	45.89
31566	31018	7.34	47.03
32265	33083	11.68	39.18
33083	33083	0.00	44.27
33878	33083	11.81	42.73

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Triglycerides
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30647-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Triglycerides; Almond oil, refined; Amygdalae oleum raffinatum; Apricot kernel oil, refined; Arachidis oleum raffinatum; Avocado oil, refined; Avocado oleum raffinatum; Borage oil, refined, stabilised; Carthami oleum raffinatum; Codliver oil (type A); Corn oil, refined; Decyl oleate; Decylis oleas; Decylium oleinicum; Evening primrose oil, refined; Flaxseed oil, native; Grape seed oil, refined; Iecoris aselli oleum A; Jojoba oil, native; Laurel oil; Lecithinum e sojae, liquidus; Lini oleum virginale; Oenotherae oleum raffinatum; Oleum amygdalarum raffinatum; Oleum arachidis raffinatum; Oleum boraginis raffinatum; Oleum carthami; Oleum cucurbita semen; Oleum helianthi raffinatum; Oleum jecoris (Typ A); Oleum juglandis raffinatum; Oleum lauri expressum; Oleum lini; Oleum maydis raffinatum; Oleum oenotherae raffinatum; Oleum olivarum nativum; Oleum omega-3 18/12; Oleum persea americana raffinatum; Oleum pruni armeniaca raffinatum; Oleum rapae raffinatum; Oleum sesami raffinatum; Oleum simmondsia chinensis; Oleum sojae raffinatum; Oleum tritici raffinatum; Oleum vitis viniferae raffinatum; Oleyl oleate (stabilised); Oleylis oleas; Oleylium oleinicum; Olivae oleum virginale; Olive oil, native; Omega-3-fish oil; Peanut oil, refined; Pumpkin seed oil, pressed; Rapae oleum raffinatum; Rapeseed oil, refined; Safflower oil, refined type I; Sesame oil, refined; Soy lecithin, liquid; Soy oil, refined; Sunflower oil, refined; Tritici aestivi oleum raffinatum; Walnut oil, refined; Wheatgerm oil, refined

Special notes

When selecting the *Triglycerides* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Triglycerides*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Rapeseed oil, re...	12054402	30647	40	20120321
Henry Lamotte	Omega-3-fish oil	0123108	30679	40	20120613
Caelo	Oleyl oleate (st...	12220205	30754	40	20120803
Caelo	Apricot kernel o...	11346402	30840	40	20111214
Caelo	Evening primrose...	12325103	30855	40	20121107
Caelo	Oleyl oleate (st...	12220207	30863	40	20120803
Caelo	Sesame oil, refi...	12054508	30873	40	20120323
Caelo	Codliver oil (ty...	12062501	30878	40	20120411
Caelo	Safflower oil, r...	13011302	30981	80	1407318
Caelo	Corn oil, refined	12141511	31016	40	1407320
Euro OTC	Soy lecithin, li...	1305019-01	31034	40	20130603
Caelo	Grape seed oil, ...	13065302	31121	40	1407373
Caelo	Laurel oil	13024511	31122	40	1403397
Caelo	Borage oil, refi...	13170606	31161	40	1407319
Bombastus	Flaxseed oil, na...	281489	31172	40	1407202
Henry Lamotte	Omega-3-fish oil	0123780	31178	40	20130403
Caelo	Wheatgerm oil, r...	13011508	31210	40	1407215
Caelo	Rapeseed oil, re...	13011008	31215	40	1406536
Caelo	Walnut oil, refi...	13202401	31216	40	1406537
Caelo	Sunflower oil, r...	13011411	31217	40	1406538
Caelo	Peanut oil, refi...	13306106	31246	40	1406534
Caelo	Almond oil, refi...	13350007	31250	40	1407348
Caelo	Corn oil, refined	13249005	31274	40	1407347
Caelo	Jojoba oil, native	13287709	31322	40	AR-14-FG-012487-01
Caelo	Flaxseed oil, na...	14087903	31326	40	20140401
Caelo	Pumpkin seed oil...	13358504	31332	40	AR-14-FG-012336-02
Caelo	Avocado oil, ref...	13361406	31336	40	AR-15-FG-002721-01
Caelo	Almond oil, refi...	14042704	31337	40	20140211
Caelo	Corn oil, refined	13249005	31344	40	20130819
Caelo	Peanut oil, refi...	14115404	31346	40	20140415
Caelo	Sesame oil, refi...	13397402	31350	40	AR-14-FG-012338-01
Caelo	Olive oil, native	14024409	31356	40	20140130
Caelo	Codliver oil (ty...	13396403	31361	40	AR-14-FG-012167-01
Caelo	Safflower oil, r...	13358705	31362	40	20131122
Caelo	Sunflower oil, r...	13225110	31374	40	20130813
Caelo	Apricot kernel o...	13357604	31375	40	1412121
Caelo	Rapeseed oil, re...	13236506	31376	40	20130806
Caelo	Soy oil, refined	13361305	31378	40	20131122
Caelo	Walnut oil, refi...	13364501	31380	40	20131122
Caelo	Decyl oleate	13378403	31382	40	AR-15-FG-000232-01
Caelo	Borage oil, refi...	14161904	31384	40	20140522
Caelo	Evening primrose...	14010915	31460	40	AR-14-FG-012898-01
Caelo	Wheatgerm oil, r...	13261604	31491	40	1412119
Caelo	Oleyl oleate (st...	14005001	31496	40	AR-15-FG-001552-01
Euro OTC	Soy lecithin, li...	1404022-01	31523	40	8779
Henry Lamotte	Omega-3-fish oil	0124495	31540	40	20140325
Caelo	Evening primrose...	14258610	31572	40	AR-15-FG-002720-01
Euro OTC	Decyl oleate	1309018-01	31584	40	20131007
Euro OTC	Olive oil, native	1407052-02	31585	40	20140812
Euro OTC	Evening primrose...	1407002-02	31589	40	20140718
Caelo	Pumpkin seed oil...	14199402	31609	40	1507428
Caelo	Avocado oil, ref...	14363703	31662	40	20150114
Fagron	Codliver oil (ty...	14D16-B08-304341	31677	40	20140520
Euro OTC	Evening primrose...	1506001	31708	40	20150616*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Henry Lamotte	Omega-3-fish oil	1000000642	31729	40	1601423
Euro OTC	Almond oil, refi...	1510014	31770	40	20151022*
Euro OTC	Evening primrose...	1509005	31775	40	20150921*
Henry Lamotte	Omega-3-fish oil	1000000642	31792	40	20150414
Euro OTC	Avocado oil, refi...	1501015	31801	40	20150121*
Fagron	Codliver oil (ty...	15A07-B05-304428	31825	40	20150116
Caelo	Laurel oil	15216401	31833	30	1606485
Caelo	Soy oil, refined	15239302	31834	40	20150817
Henry Lamotte	Omega-3-fish oil	1000001175	31875	40	20150626
Caelo	Olive oil, native	161446	32195	40	20160524*
Caelo	Peanut oil, refi...	162495	32680	40	20160831*
Caelo	Apricot kernel o...	162361	32683	40	20160818*
Caelo	Wheatgerm oil, r...	16137403	32730	40	20160708
Caelo	Sunflower oil, r...	15422804	32745	40	20151218
Caelo	Grape seed oil, ...	16092402	32989	40	20160407
Bombastus	Olive oil, native	306949	32991	40	20170116
Caelo	Peanut oil, refi...	170428	33059	40	20170215*
Caelo	Apricot kernel o...	170555	33060	30	20170308*
Caelo	Flaxseed oil, na...	170467	33064	40	20170222*
Caelo	Avocado oil, refi...	170388	33067	40	20170209*
Caelo	Decyl oleate	170595	33070	40	20170308*
Caelo	Jobaba oil, native	16070602	33254	40	20160316
Caelo	Laurel oil	16275701	33255	40	20161005
Caelo	Soy oil, refined	171217	33309	40	20170509*
Caelo	Sesame oil, refi...	171777	33313	40	20170620*
Caelo	Borage oil, refi...	16266007	33510	40	20161027
Henry Lamotte	Jobaba oil, native	8002391001	33876	40	20161124
Caelo	Walnut oil, refi...	173248	33916	40	20170512*
Caelo	Pumpkin seed oil...	171379	33917	40	20170606*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 3340 spectra of 83 reference samples from the substance/substance group *Triglycerides*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 81 different batches.
- 32 204 spectra from a total of 759 batches from further 187 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 4921 spectra from 129 reference samples of the substance/substance group *Triglycerides*.
- These include spectra of independent samples from 95 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Olive oil, native	90417099	30044	9
Caelo	Peanut oil, refined	90753159	30049	10
Caelo	Borage oil, refined, stabili...	09267308	30162	40
Caelo	Laurel oil	10185405	30168	40
Caelo	Grape seed oil, refined	10066303	30192	29
Caelo	Sunflower oil, refined	10159002	30200	37
Bombastus	Olive oil, native	00254963	30226	20
Caelo	Peanut oil, refined	1081502	30227	10
Caelo	Olive oil, native	10314504	30249	40
Caelo	Almond oil, refined	11032604	30257	40
Caelo	Almond oil, refined	11032604	30259	40
Caelo	Avocado oil, refined	10215207	30270	40
Caelo	Corn oil, refined	10104608	30271	40
Caelo	Peanut oil, refined	11000503	30272	130
Caelo	Olive oil, native	10314508	30273	130
Caelo	Almond oil, refined	11032605	30274	130
Caelo	Avocado oil, refined	10215207	30275	130
Caelo	Olive oil, native	11063901	30288	50
Caelo	Avocado oil, refined	10215207	30290	32
Bombastus	Jobba oil, native	00256221	30294	50
Caelo	Jobba oil, native	00252449	30314	40
Euro OTC	Soy lecithin, liquid	1006017-01	30321	30
Caelo	Peanut oil, refined	11163201	30331	40
Caelo	Avocado oil, refined	00243981	30347	40
Henry Lamotte	Peanut oil, refined	5052601	30366	30
Statfold	Flaxseed oil, native	17732	30413	40
Caelo	Walnut oil, refined	10328802	30419	40
Bombastus	Flaxseed oil, native	00269324	30522	40
Caelo	Wheatgerm oil, refined	11227305	30550	40
Caelo	Peanut oil, refined	14115404	31346 [†]	20
Caelo	Sesame oil, refined	13397402	31350 [†]	20
Caelo	Olive oil, native	14024409	31356 [†]	20
Caelo	Codliver oil (type A)	13396403	31361 [†]	20
Caelo	Safflower oil, refined type I	13358705	31362 [†]	20
Caelo	Evening primrose oil, refined	14010915	31460 [†]	20
Caelo	Jobba oil, native	13287709	31482	60
Caelo	Oleyl oleate (stabilised)	14005001	31496 [†]	20
Caelo	Jobba oil, native	13287709	31519	60
Euro OTC	Soy lecithin, liquid	1404022-01	31523 [†]	20
Henry Lamotte	Omega-3-fish oil	0124495	31540 [†]	20
Caelo	Evening primrose oil, refined	14258610	31572 [†]	20

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Decyl oleate	1309018-01	31584 [†]	20
Euro OTC	Olive oil, native	1407052-02	31585 [†]	20
Euro OTC	Evening primrose oil, refined	1407002-02	31589 [†]	20
Caelo	Pumpkin seed oil, pressed	14199402	31609 [†]	20
Henry Lamotte	Codliver oil (type A)	0125348	31643	60
Caelo	Laurel oil	14176301	31657	15
Caelo	Avocado oil, refined	14363703	31662 [†]	20
Fagron	Codliver oil (type A)	14D16-B08-304341	31677 [†]	20
Euro OTC	Evening primrose oil, refined	1506001	31708 [†]	20
Henry Lamotte	Omega-3-fish oil	1000000642	31729 [†]	20
Caelo	Laurel oil	15216402	31745	60
Euro OTC	Almond oil, refined	1510014	31770 [†]	20
Euro OTC	Evening primrose oil, refined	1509005	31775 [†]	20
Henry Lamotte	Omega-3-fish oil	1000000642	31792 [†]	20
Euro OTC	Avocado oil, refined	1501015	31801 [†]	20
Fagron	Codliver oil (type A)	15A07-B05-304428	31825 [†]	20
Caelo	Laurel oil	15216401	31833 [†]	30
Caelo	Soy oil, refined	15239302	31834 [†]	20
Henry Lamotte	Omega-3-fish oil	1000001175	31875 [†]	20
Euro OTC	Olive oil, native	1602022	32186	40
Euro OTC	Almond oil, refined	1603025	32188	40
Caelo	Decyl oleate	161071	32194	40
Caelo	Avocado oil, refined	161047	32197	40
Caelo	Almond oil, refined	161444	32198	40
Caelo	Walnut oil, refined	161007	32199	40
Caelo	Oleyl oleate (stabilised)	161049	32202	40
Caelo	Peanut oil, refined	161195	32206	40
Caelo	Avocado oil, refined	161859	32315	40
Caelo	Rapeseed oil, refined	161850	32316	40
Caelo	Safflower oil, refined type I	161849	32334	40
Euro OTC	Avocado oil, refined	1605034	32457	40
Euro OTC	Evening primrose oil, refined	1605027	32463	40
Caelo	Evening primrose oil, refined	162813	32682	40
Caelo	Borage oil, refined, stabili...	162660	32685	40
Caelo	Corn oil, refined	162603	32689	40
Caelo	Codliver oil (type A)	162739	32690	40
Caelo	Sunflower oil, refined	162758	32691	40
Euro OTC	Olive oil, native	1608006	32709	40
Caelo	Soy oil, refined	15239302	32729	40
Henry Lamotte	Omega-3-fish oil	1000002824	32734	40
Caelo	Grape seed oil, refined	16092401	32746	40
Caelo	Pumpkin seed oil, pressed	15410004	32783	40
Caelo	Apricot kernel oil, refined	16236106	32931	40
Bombastus	Olive oil, native	306949	32990	40
Caelo	Olive oil, native	163274	33063	30
Caelo	Almond oil, refined	170422	33066	40
Caelo	Olive oil, native	170387	33068	40
Caelo	Corn oil, refined	170550	33069	40
Caelo	Sesame oil, refined	163125	33075	40
Caelo	Peanut oil, refined	163553	33085	40
Euro OTC	Evening primrose oil, refined	1703020	33227	40
Euro OTC	Almond oil, refined	1701036	33231	40
Caelo	Olive oil, native	171216	33306	40
Caelo	Codliver oil (type A)	171378	33308	40
Caelo	Almond oil, refined	171765	33314	40
Caelo	Peanut oil, refined	171690	33315	40
Caelo	Soy oil, refined	16347104	33462	39
Henry Lamotte	Sunflower oil, refined	8002491005	33472	40
Caelo	Grape seed oil, refined	16339703	33475	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Grape seed oil, refined	16339703	33476	40
Caelo	Walnut oil, refined	16100703	33488	40
Caelo	Wheatgerm oil, refined	16321301	33489	40
Caelo	Apricot kernel oil, refined	16236109	33512	40
Caelo	Sesame oil, refined	16312505	33618	40
Caelo	Safflower oil, refined type I	17092201	33653	40
Caelo	Decyl oleate	172156	33671	40
Caelo	Rapeseed oil, refined	172139	33673	40
Caelo	Olive oil, native	172365	33674	40
Caelo	Flaxseed oil, native	172364	33677	40
Caelo	Apricot kernel oil, refined	172518	33700	40
Euro OTC	Soy lecithin, liquid	1703001-01	33838	40
Euro OTC	Evening primrose oil, refined	1710011	33849	40
Caelo	Borage oil, refined, stabili...	17055407	33853	40
Bombastus	Flaxseed oil, native	309176	33880	40
Caelo	Jojoba oil, native	170733	33892	40
Caelo	Wheatgerm oil, refined	172903	33915	40
Caelo	Oleyl oleate (stabilised)	173581	33920	40
Euro OTC	Almond oil, refined	1801030	33924	40
Caelo	Laurel oil	18028202	33948	40
Euro OTC	Olive oil, native	1802011	34005	40
Caelo	Pumpkin seed oil, pressed	180837	34061	40
Caelo	Safflower oil, refined type I	180999	34062	40
Euro OTC	Soy lecithin, liquid	1803030	34094	40
Euro OTC	Evening primrose oil, refined	1803006	34095	40
Henry Lamotte	Omega-3-fish oil	1000007067	34124	40
Euro OTC	Olive oil, native	1805023	34183	40
Euro OTC	Almond oil, refined	1807007	34195	40
Euro OTC	Olive oil, native	1808026	34215	40

- 88 953 spectra from a total of 2045 batches from further 412 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 8294 spectra from 453 *Apo-Ident* customers from 2044 batches from the substance/substance group *Triglycerides*.
- These include spectra of independent samples from 1986 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
VDL 21.05;Preis 9,84	Almond oil, refined	13236406//05R05	1
Caelo	Almond oil, refined	13350003	14

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Almond oil, refined	13350001	18
Caelo	Almond oil, refined	14139305	2
Caelo	Almond oil, refined	14042703	7
Caelo	Almond oil, refined	14139310	14
Caelo	Almond oil, refined	14139311	9
Caelo	Almond oil, refined	14215806	10
Caelo	Almond oil, refined	14215816	1
Caelo	Almond oil, refined	14215807	2
Caelo	Almond oil, refined	15078404	11
Caelo	Almond oil, refined	15409712	6
Bombastus	Almond oil, refined	283800	2
Bombastus	Almond oil, refined	282093	3
Bombastus	Almond oil, refined	284967	3
Bombastus	Almond oil, refined	286253	2
Bombastus	Almond oil, refined	288471	3
Bombastus	Almond oil, refined	287416	1
Bombastus	Almond oil, refined	295615	2
Lamotte/Gehe	Almond oil, refined	5615001	1
Henry Lamotte/Sanacorp	Almond oil, refined	130145	1
Caelo	Almond oil, refined	12153502	6
Lamotte / Gehe	Almond oil, refined	5748601	4
Sanacorp	Almond oil, refined	12357103	1
Caelo	Almond oil, refined	12357111	12
Caelo	Almond oil, refined	12357112	3
Caelo	Almond oil, refined	13178002	14
Caelo	Almond oil, refined	13178004	15
Sanacorp	Almond oil, refined	13178004	1
Caelo	Almond oil, refined	13236402	19
Bombastus	Almond oil, refined	277971	5
Bombastus	Almond oil, refined	281040	5
Fagron	Almond oil, refined	281040	2
Henry Lamotte Oils	Almond oil, refined	5438302	5
Lamotte/Sanacorp	Almond oil, refined	5438302	2
Henry Lamotte/Gehe	Almond oil, refined	5415901	2
Fagron	Almond oil, refined	5438301	2
Henry Lamotte/N/PZN:731...	Almond oil, refined	5438302	2
Caelo	Apricot kernel oil, refined	13253203	2
Primavera	Apricot kernel oil, refined	3612131	1
Caelo	Apricot kernel oil, refined	13011804	1
taosis/Phönix	Apricot kernel oil, refined	134933-124392	1
Caelo	Apricot kernel oil, refined	13011805	1
Bombastus	Apricot kernel oil, refined	274930	1
Caelo	Peanut oil, refined	13358803	15
Fagron	Peanut oil, refined	13G30-N05	8
Caelo	Peanut oil, refined	14115406	17
Caelo	Peanut oil, refined	14115405	14
Caelo	Peanut oil, refined	14246901	16
Caelo	Peanut oil, refined	14355904	14
Caelo	Peanut oil, refined	12203803	19
Caelo	Peanut oil, refined	15015103	14
Euro OTC	Peanut oil, refined	5592601	1
Caelo	Peanut oil, refined	12274302	29
Caelo	Peanut oil, refined	12274305	3
Caelo	Peanut oil, refined	12323205	3
Sanacorp	Peanut oil, refined	13129203	1
Caelo	Peanut oil, refined	13306103	27
Caelo	Peanut oil, refined	13306102	8
Bombastus	Peanut oil, refined	284286	5
Caelo	Avocado oil, refined	12248312	3
Euro OTC	Avocado oil, refined	1405007-01	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
VDL 06.05.Preis 6,70	Avocado oil, refined	700325//05R06	1
Caelo	Avocado oil, refined	13129108	7
Caelo	Avocado oil, refined	13361407	2
Euro OTC	Avocado oil, refined	1301017-01	1
Sanacorp	Avocado oil, refined	1304039-01	1
Euro OTC	Avocado oil, refined	1304039-01	6
Euro OTC	Avocado oil, refined	1305046-01	5
Caelo	Avocado oil, refined	13129103	7
Henry Lamotte Oils GmbH	Avocado oil, refined	1000001098	1
Caelo	Avocado oil, refined	15155206	7
Henry Lamotte/AHD	Avocado oil, refined	1000004363	1
Lamotte GmbH	Avocado oil, refined	1000001451	1
Caelo	Avocado oil, refined	15268203	10
Caesar & Loretz/Noweda	Avocado oil, refined	15268203	1
Caelo	Avocado oil, refined	15268205	5
Caelo	Avocado oil, refined	15268213	2
Caelo	Avocado oil, refined	15426303	2
Caelo	Avocado oil, refined	15426304	5
Caelo	Avocado oil, refined	15426305	8
Caelo	Avocado oil, refined	15426306	2
Euro OTC	Avocado oil, refined	1605034-01	6
Euro OTC	Avocado oil, refined	1605034-03	1
Caelo	Avocado oil, refined	16104702	14
Caelo	Avocado oil, refined	16104703	3
Caelo	Avocado oil, refined	16185907	5
Caelo	Avocado oil, refined	16185908	4
Caelo	Avocado oil, refined	17038802	9
Caelo	Avocado oil, refined	1618908	1
phönix	Avocado oil, refined	17038802	1
Caesar & Loretz/Noweda	Avocado oil, refined	17038802	1
Wepa/ Phoenix	Avocado oil, refined	17038802	1
Caesar & Loretz/Noweda	Avocado oil, refined	17038803	1
Caelo	Avocado oil, refined	18044901	1
Caelo	Avocado oil, refined	18007501	1
Bombastus	Avocado oil, refined	294183	4
Caelo	Avocado oil, refined	2515I-01960	1
Caelo	Avocado oil, refined	300403	1
Ebert u. Jacobi 16.0...	Avocado oil, refined	302706	1
Bombastus	Avocado oil, refined	300403	1
Bombastus	Avocado oil, refined	306789	1
Spinnrad GmbH	Avocado oil, refined	916301202	1
Henry Lamotte Oils GmbH	Avocado oil, refined	125529	2
Lamotte/Sanacorp	Avocado oil, refined	130214	1
Henry lamotte Oils GmbH...	Avocado oil, refined	130243	1
Henry Lamotte Oils GmbH	Avocado oil, refined	130243	1
Euro OTC	Avocado oil, refined	2041404	1
Henry Lamotte/Sana	Avocado oil, refined	1000001098	1
Henry Lamotte/Anzag	Avocado oil, refined	1000001098	2
Euro OTC	Avocado oil, refined	13111404	1
Phoenix 19.03.2014	Avocado oil, refined	13129102	1
Nachprüfung	Avocado oil, refined	13129103	2
Caelo	Avocado oil, refined	13361406	6
Caelo	Avocado oil, refined	13129111	1
Phoenix 29.04.2014	Avocado oil, refined	13361406	1
Holdermann	Avocado oil, refined	13361406	1
Großhandel	Avocado oil, refined	13361407	1
Caelo	Avocado oil, refined	13361411	2
Caelo	Avocado oil, refined	13361414	7
Caesar & Loretz GmbH Hi...	Avocado oil, refined	13361414	1
Audor Pharma	Avocado oil, refined	13361414	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Avocado oil, refined	13361441	1
Fagron	Avocado oil, refined	13J25-B03	1
Raben-Apotheke Mildenau	Avocado oil, refined	1402027-01	1
Euro OTC	Avocado oil, refined	1402027-01	7
Bombastus	Avocado oil, refined	1402027-01	1
Euro OTC	Avocado oil, refined	1410027-01	4
Sanacorp	Avocado oil, refined	14170104	1
Caelo	Avocado oil, refined	14170104	11
Caelo	Avocado oil, refined	14170105	1
Caelo	Avocado oil, refined	14363702	10
Nachprüfung	Avocado oil, refined	143637021	1
Caelo	Avocado oil, refined	14363703	9
Caelo	Avocado oil, refined	14363704	1
Caelo	Avocado oil, refined	14363705	4
Caelo	Avocado oil, refined	14363706	2
Caelo	Avocado oil, refined	1436705	1
Euro OTC	Avocado oil, refined	1501015-01	3
Euro OTC	Avocado oil, refined	1501015V-01	1
Caelo	Avocado oil, refined	15155203	2
Caelo	Avocado oil, refined	15155207	2
phönix	Avocado oil, refined	16185908	1
Caelo	Avocado oil, refined	19081501	1
Euro OTC	Avocado oil, refined	28041506	1
Bombastus	Avocado oil, refined	286519	1
Bombastus	Avocado oil, refined	289240	1
Euro OTC	Avocado oil, refined	30121424	1
Bombastus	Avocado oil, refined	267660	1
Caelo	Avocado oil, refined	3413M-01960	2
Bombastus	Avocado oil, refined	268579	1
Henry Lamotte Oils	Avocado oil, refined	122470	1
Henry Lamotte Oils GmbH	Avocado oil, refined	122490	2
Lamotte	Avocado oil, refined	122490	1
Henry Lamotte	Avocado oil, refined	130143	1
HenryLamotte/Phoenix	Avocado oil, refined	130143	1
Henry Lamotte Oils GmbH	Avocado oil, refined	130143	2
Lamotte	Avocado oil, refined	130143	1
Caelo	Avocado oil, refined	130172	1
Lamotte/Sanacorp	Avocado oil, refined	130193	1
Caelo	Avocado oil, refined	2081303	1
Bombastus	Avocado oil, refined	2091235000130	1
PVL	Avocado oil, refined	2410131	1
Euro OTC	Avocado oil, refined	1112018-01	1
PVL	Avocado oil, refined	2807133	1
Caelo	Avocado oil, refined	11147607	1
Caelo	Avocado oil, refined	11227207	2
Anzag 16.10.12 ek: 3,90EUR	Avocado oil, refined	1204047-01	1
Caesar/Noweda	Avocado oil, refined	12064801	1
Großhandel	Avocado oil, refined	12064801	1
Caelo	Avocado oil, refined	12064802	2
Caelo	Avocado oil, refined	12064804	2
Caelo	Avocado oil, refined	12064810	3
Caesar/Noweda	Avocado oil, refined	12064812	2
Caelo	Avocado oil, refined	12064813	3
Anzag 06.02.2013	Avocado oil, refined	1208028-01	1
Caelo	Avocado oil, refined	12248302	2
Großhandel	Avocado oil, refined	12248302	1
Caelo	Avocado oil, refined	12248310	11
Caelo	Avocado oil, refined	12248303	2
Phoenix 30.04.2013	Avocado oil, refined	12248310	2
Caesar & Loretz/ Noweda	Avocado oil, refined	12248312	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Avocado oil, refined	12248315	2
Caesar & Loretz GmbH	Avocado oil, refined	12H24-N02	1
Euro OTC	Avocado oil, refined	1301039-01	1
Fagron	Avocado oil, refined	12H24-N02	1
Caelo	Avocado oil, refined	13031303	1
Caelo	Avocado oil, refined	1304039-1	1
gehe 15.01.2014	Avocado oil, refined	1305046-01	1
Großhandel	Avocado oil, refined	13129102	1
Noveda	Avocado oil, refined	13129103	1
Caelo	Avocado oil, refined	13129104	2
Caelo	Avocado oil, refined	13129106	3
spinnrad	Avocado oil, refined	212301211	1
Caelo	Avocado oil, refined	28111308	1
Bombastus	Avocado oil, refined	274935	1
Phönix 15.01.2014	Avocado oil, refined	279025	1
EuRho	Avocado oil, refined	4111M-01960	1
spinnrad	Avocado oil, refined	413301212	1
Caelo	Avocado oil, refined	4718201	2
Caelo	Avocado oil, refined	5216901	1
Taoasis	Avocado oil, refined		1
Taoasis	Avocado oil, refined	92010-105977	1
Caelo	Almond oil, refined	2204130	1
Lamotte/Alliance	Borage oil, refined, stabili...	125354	1
Henry Lamotte	Borage oil, refined, stabili...	145579	1
Henry Lamotte Oils GmbH	Borage oil, refined, stabili...	130402	1
Henry Lamotte/Anzag	Borage oil, refined, stabili...	125354	1
Henry Lamotte/Anzag	Borage oil, refined, stabili...	1000001388	1
Caelo	Borage oil, refined, stabili...	13170607	1
Caelo	Borage oil, refined, stabili...	13372104	6
Spinnrad	Borage oil, refined, stabili...	1115204582	1
Caelo	Borage oil, refined, stabili...	13372106	2
Caelo	Borage oil, refined, stabili...	13372107	8
Caelo	Borage oil, refined, stabili...	14161901	4
Caelo	Borage oil, refined, stabili...	14161904	4
Caelo	Borage oil, refined, stabili...	14161906	2
Caelo	Borage oil, refined, stabili...	1416906	1
Caelo	Borage oil, refined, stabili...	14203001	3
Caelo	Borage oil, refined, stabili...	14203004	11
Phönix 14.01.15	Borage oil, refined, stabili...	14203004	1
Caelo	Borage oil, refined, stabili...	14203005	12
Sanacorp	Borage oil, refined, stabili...	14203006	1
Caelo	Borage oil, refined, stabili...	14203009	1
Caelo	Borage oil, refined, stabili...	14203011	5
Caelo	Borage oil, refined, stabili...	14203012	1
Euro OTC	Borage oil, refined, stabili...	14203012	1
Caelo	Borage oil, refined, stabili...	1461906	1
Caelo	Borage oil, refined, stabili...	15063401	2
Fagron	Borage oil, refined, stabili...	14A20-B03	1
Caelo	Borage oil, refined, stabili...	15063402	12
Caelo	Borage oil, refined, stabili...	15063406	2
Caelo	Borage oil, refined, stabili...	15063407	14
Euro OTC	Borage oil, refined, stabili...	15063407	1
Caelo	Borage oil, refined, stabili...	15063408	3
Caelo	Borage oil, refined, stabili...	15063409	2
Caelo	Borage oil, refined, stabili...	15063411	5
Caesar & Loretz GmbH / ...	Borage oil, refined, stabili...	15063412	1
Caelo	Borage oil, refined, stabili...	15063412	2
Caelo	Borage oil, refined, stabili...	15063413	2
Caelo	Borage oil, refined, stabili...	15305001	12
Caelo	Borage oil, refined, stabili...	15305003	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Borage oil, refined, stabili...	15305005	3
Caelo	Borage oil, refined, stabili...	16266007	10
Caelo	Borage oil, refined, stabili...	26041404	1
PG-Naturpharm	Borage oil, refined, stabili...	397-045	2
Aurica Vital	Borage oil, refined, stabili...	28/05/15	1
Caelo	Borage oil, refined, stabili...	6170717	1
Lamotte/Noweda	Borage oil, refined, stabili...	5404602	1
Aurica/Noweda	Borage oil, refined, stabili...	957-038	1
Caelo	Borage oil, refined, stabili...	722915305009	1
Caelo	Borage oil, refined, stabili...	1000006595	1
Nachprüfung	Borage oil, refined, stabili...	1503412	1
Caelo	Borage oil, refined, stabili...	15305002	2
Stifts-Apotheke	Borage oil, refined, stabili...	15305002	1
Caelo	Borage oil, refined, stabili...	15305009	7
Caelo	Borage oil, refined, stabili...	15305011	1
Caelo	Borage oil, refined, stabili...	1534-08	1
Caelo	Borage oil, refined, stabili...	16052301	7
Caelo	Borage oil, refined, stabili...	16052302	2
Caelo	Borage oil, refined, stabili...	16052307	6
Caelo	Borage oil, refined, stabili...	16266002	4
Caelo	Borage oil, refined, stabili...	16266011	5
Caelo	Borage oil, refined, stabili...	17055401	5
Sanacorp W.E: 01.02.18 ...	Borage oil, refined, stabili...	17307801	2
Caelo	Borage oil, refined, stabili...	17307801	6
Euro OTC	Borage oil, refined, stabili...	17307801	1
Caelo	Borage oil, refined, stabili...	17361804	3
Caelo	Borage oil, refined, stabili...	18080501	2
Henry Lamotte GmbH	Borage oil, refined, stabili...	5594400001	1
Spinnrad GmbH	Borage oil, refined, stabili...	616204581	1
Henry Lamotte/Anzag	Borage oil, refined, stabili...	5594400001	1
Aurica Vital	Borage oil, refined, stabili...	Z13902/A	1
Aurica	Borage oil, refined, stabili...	Z13902/A	1
Caelo	Borage oil, refined, stabili...	11260102	2
Caelo	Borage oil, refined, stabili...	123118	1
	Borage oil, refined, stabili...	123387	1
Henry Lamotte GmbH Bremen	Borage oil, refined, stabili...	122791	1
Caelo	Borage oil, refined, stabili...	11002607	1
Caelo	Borage oil, refined, stabili...	11260106	2
Caelo	Borage oil, refined, stabili...	12027004	1
Gehe	Borage oil, refined, stabili...	12027008	1
Caelo	Borage oil, refined, stabili...	12027008	2
Noweda	Borage oil, refined, stabili...	12161202	2
Caelo	Borage oil, refined, stabili...	12152108	1
Noweda	Borage oil, refined, stabili...	12027006	1
Caelo	Borage oil, refined, stabili...	12161202	8
Kehr	Borage oil, refined, stabili...	12161202	1
Caesar/Noweda	Borage oil, refined, stabili...	12161202	1
Gehe	Borage oil, refined, stabili...	12161203	1
Caelo	Borage oil, refined, stabili...	12161205	3
Caelo	Borage oil, refined, stabili...	12161203	1
Caelo	Borage oil, refined, stabili...	12161207	5
Noweda	Borage oil, refined, stabili...	12355203	1
Caelo	Borage oil, refined, stabili...	1235504	1
Caelo	Borage oil, refined, stabili...	12161209	1
Caelo	Borage oil, refined, stabili...	12355203	10
Phoenix	Borage oil, refined, stabili...	12355203	1
Caelo	Borage oil, refined, stabili...	12355205	4
Caelo	Borage oil, refined, stabili...	1235523	1
Caelo	Borage oil, refined, stabili...	12365203	1
J.H.Müller/direkt/17.05...	Borage oil, refined, stabili...	130302	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar/Noweda	Borage oil, refined, stabili...	13136501	1
Noweda	Borage oil, refined, stabili...	13136504	1
Caelo	Borage oil, refined, stabili...	13136504	10
Caesar/Noweda	Borage oil, refined, stabili...	13136504	1
Caelo	Borage oil, refined, stabili...	1316504	1
Caelo	Borage oil, refined, stabili...	13170601	4
Caelo	Borage oil, refined, stabili...	13170605	3
Caelo	Borage oil, refined, stabili...	13170604	1
Caelo	Borage oil, refined, stabili...	13170606	18
Noweda	Borage oil, refined, stabili...	13170606	1
	Borage oil, refined, stabili...	13170606	1
Caelo	Borage oil, refined, stabili...	13170610	7
Caelo	Borage oil, refined, stabili...	13170613	4
Caesar/Noweda	Borage oil, refined, stabili...	13372103	1
Noweda	Borage oil, refined, stabili...	13372101	1
Noweda	Borage oil, refined, stabili...	13372104	1
Phoenix	Borage oil, refined, stabili...	14161906	1
Caelo	Borage oil, refined, stabili...	1713M-02881	1
Henry Lamotte	Borage oil, refined, stabili...	5402901	3
Anzag 15.10.12	Borage oil, refined, stabili...	965-028	2
Bergland-Pharma	Borage oil, refined, stabili...	H103019	2
Henry Lamotte Oils/Sana...	Borage oil, refined, stabili...	5402901	1
Caelo	Oleyl oleate (stabilised)	15321701	9
Caelo	Oleyl oleate (stabilised)	17013008	1
Noweda	Oleyl oleate (stabilised)	15321706	1
Fagron	Oleyl oleate (stabilised)	7858448	1
Caelo	Oleyl oleate (stabilised)	25062015A	1
Caelo	Oleyl oleate (stabilised)	11292610	1
Caelo	Oleyl oleate (stabilised)	10298401	2
Fagron	Oleyl oleate (stabilised)	11108-N03	1
Caelo	Oleyl oleate (stabilised)	12067701	1
Caelo	Oleyl oleate (stabilised)	12067704	9
Caelo	Oleyl oleate (stabilised)	12067705	2
Sanacorp	Oleyl oleate (stabilised)	12067705	1
Caelo	Oleyl oleate (stabilised)	12220201	1
Caelo	Oleyl oleate (stabilised)	12142803	1
VDL	Oleyl oleate (stabilised)	12220205	2
Caelo	Oleyl oleate (stabilised)	12220204	1
Caelo	Oleyl oleate (stabilised)	12220205	5
Caelo	Oleyl oleate (stabilised)	12220207	3
Gehe	Oleyl oleate (stabilised)	12220207	1
Caelo	Oleyl oleate (stabilised)	12220210	1
Caelo	Oleyl oleate (stabilised)	13167109	4
Caelo	Oleyl oleate (stabilised)	13096407	1
Caelo	Oleyl oleate (stabilised)	13167107	1
Caelo	Oleyl oleate (stabilised)	13144401	1
Caelo	Oleyl oleate (stabilised)	13096401	1
Caelo	Oleyl oleate (stabilised)	122220210	1
gehe 20.08.2013	Oleyl oleate (stabilised)	13167109	1
Caelo	Oleyl oleate (stabilised)	13291607	1
Caelo	Oleyl oleate (stabilised)	14005001	4
Caelo	Oleyl oleate (stabilised)	21061302	1
SysKem Chemie	Oleyl oleate (stabilised)	35885	1
Caelo	Oleyl oleate (stabilised)	26071210	1
SysKem Chemie	Oleyl oleate (stabilised)	8002060413	1
Caelo	Oleyl oleate (stabilised)		2
Caelo	Decyl oleate	15082206	3
Caelo	Decyl oleate	15191502	3
Caelo	Decyl oleate	15248405	7
Caelo	Decyl oleate	16107101	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Decyl oleate	16224309	1
Caelo	Decyl oleate	13378403	2
Caelo	Decyl oleate	14005001	2
Caelo	Decyl oleate	14005002	1
Sanacorp	Decyl oleate	13378403	1
Raben-Apotheke Mildenau	Decyl oleate	14814-B03-292195	1
Caelo	Decyl oleate	14254603	1
Caelo	Olive oil, native	13350103	13
Caelo	Decyl oleate	15082205	1
Caelo	Decyl oleate	15191501	1
Sanacorp	Decyl oleate	15082206	1
Caelo	Decyl oleate	13291606	2
Euro OTC	Decyl oleate	1309018-01	1
Caelo	Decyl oleate	15047702	2
Caelo	Decyl oleate	15248408	1
Caelo	Decyl oleate	15393704	1
Caelo	Decyl oleate	16223409	1
Caelo	Decyl oleate	16224302	2
Fagron	Peanut oil, refined	133498	1
Fagron	Peanut oil, refined	1410-B07-300438	1
Fagron	Peanut oil, refined	14110-B07-300438	2
Caelo	Peanut oil, refined	15273402	14
Nachprüfung	Peanut oil, refined	1426903	1
Caelo	Peanut oil, refined	15273403	11
Caelo	Peanut oil, refined	15331501	7
Caelo	Peanut oil, refined	1416A-01954	1
Caelo	Peanut oil, refined	15331502	20
Caelo	Peanut oil, refined	15331503	10
Caelo	Peanut oil, refined	15331506	1
Fagron	Peanut oil, refined	15B11B02304989	2
Fagron	Peanut oil, refined	15E07-B04-308779	2
Caelo	Peanut oil, refined	16071101	11
Caelo	Peanut oil, refined	16071102	20
Caelo	Peanut oil, refined	16071103	9
Caelo	Peanut oil, refined	16071104	11
Caesar & Loretz/AHC	Peanut oil, refined	16071104	1
Jenne	Peanut oil, refined	16071104	1
Bombastus	Peanut oil, refined	16071104	1
Caelo	Peanut oil, refined	16071105	15
Caelo	Peanut oil, refined	16071106	2
Henry Lamotte Oils GmbH	Peanut oil, refined	16071105	1
Caelo	Peanut oil, refined	160711104	1
Caelo	Peanut oil, refined	16249502	25
Caesar & Loretz GmbH	Peanut oil, refined	16249502	1
Caelo	Peanut oil, refined	16259502	1
Caelo	Peanut oil, refined	16348001	17
Caelo	Peanut oil, refined	170420802	1
Caelo	Peanut oil, refined	17042802	9
Fagron	Peanut oil, refined	16A27-B04-322519	2
Caelo	Peanut oil, refined	17169001	4
Caelo	Peanut oil, refined	17169010	2
Caelo	Peanut oil, refined	17169018	3
Caelo/Noweda	Peanut oil, refined	17169010	1
Caelo	Peanut oil, refined	17169020	1
Bombastus	Peanut oil, refined	29527	1
Fagron	Peanut oil, refined	17A10-B01-332760	1
Caelo	Peanut oil, refined	17328601	6
Caesar & Loretz/ Caesar...	Peanut oil, refined	17328601	1
Bombastus	Peanut oil, refined	295272	10
Bombastus	Peanut oil, refined	299136	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Peanut oil, refined	30100	1
Bombastus	Peanut oil, refined	302017	4
Caelo	Peanut oil, refined	302017	2
Bombastus	Peanut oil, refined	303460	5
Henry Lamotte Oils GmbH	Peanut oil, refined	8000666005	3
Bombastus	Peanut oil, refined	307541	1
Henry Lamotte	Peanut oil, refined	8000666005	3
Henry Lamotte Oils GmbH	Peanut oil, refined	5904100001	2
Ebert KKH 18.10.16	Peanut oil, refined	8001417003	2
Henry Lamotte Oils GmbH	Peanut oil, refined	8001959001	1
Caelo	Peanut oil, refined	8001417003	3
Lamotte/Phönix	Peanut oil, refined	8003146001	1
Caelo	Peanut oil, refined	13051403	1
Caelo	Peanut oil, refined		3
Fagron	Peanut oil, refined	130723	1
Fagron	Peanut oil, refined	131943	1
Caelo	Peanut oil, refined	13173003	1
Fagron	Peanut oil, refined	13306103	3
Caelo	Peanut oil, refined	13358801	10
Caelo	Peanut oil, refined	13358804	21
Fagron	Peanut oil, refined	13j25b06291837	1
Fagron	Peanut oil, refined	13J25-B06-291838	4
Caelo	Peanut oil, refined	13J25-B06-291838	1
Caleo	Peanut oil, refined	14/10-B07	1
Fagron	Peanut oil, refined	14/10-B07-300426	1
Caelo	Peanut oil, refined	14033601	18
Fagron	Peanut oil, refined	14033601	1
Caelo	Peanut oil, refined	14115403	5
Noweda	Peanut oil, refined	1410B07300438	1
Caelo	Peanut oil, refined	14115404	12
Cealo/Phoenix	Peanut oil, refined	14115405	1
Caelo	Peanut oil, refined	14246903	20
Fagron	Peanut oil, refined	14246901	1
Caelo	Peanut oil, refined	14246902	5
Caelo	Peanut oil, refined	1424693	1
Hedinger	Peanut oil, refined	14355903	1
Caelo	Peanut oil, refined	14355903	6
Fagron	Peanut oil, refined	14355904	2
Caelo	Peanut oil, refined	14355905	13
Fagron	Peanut oil, refined	14D11-B05-294672	2
Fagron	Peanut oil, refined	14D11-B05-294677	4
Anzag	Peanut oil, refined	15015103	1
Fagron	Peanut oil, refined	14I-10B07-300438	1
Caelo	Peanut oil, refined	15015104	14
Caelo	Peanut oil, refined	15015106	4
caesar&loretz,phönix	Peanut oil, refined	1508506	1
Caelo	Peanut oil, refined	15208501	3
Caelo	Peanut oil, refined	15208505	15
Caelo	Peanut oil, refined	15208506	19
Phönix Gotha	Peanut oil, refined	15273403	1
Caelo	Peanut oil, refined	1531502	1
Caesar & Loretz GmbH Hi...	Peanut oil, refined	15331502	1
Fagron	Peanut oil, refined	15E07-B04	1
Cealo/AHD	Peanut oil, refined	16071101	1
caesar &loretz GmbH,Phö...	Peanut oil, refined	16071101	1
Caesar&Loretz GmbH,PHön...	Peanut oil, refined	16348001	1
Caelo	Peanut oil, refined	20021502	1
Caelo	Peanut oil, refined	23051406	1
Bombastus	Peanut oil, refined	286393	3
Bombastus	Peanut oil, refined	283866	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Peanut oil, refined	286393	2
Caelo	Peanut oil, refined	25091405	1
Bombastus	Peanut oil, refined	288510	2
Bombastus	Peanut oil, refined	291671	5
gehe	Peanut oil, refined	310007	1
Caelo	Peanut oil, refined	30121415	1
Caelo	Peanut oil, refined	3614M-01954	1
Henry Lamotte	Peanut oil, refined	5204402	1
Henry Lamotte Oils GmbH	Peanut oil, refined	5473601	1
Lamotte/Phoenix	Peanut oil, refined	5592601	1
lamotte/Gehe	Peanut oil, refined	5473601	2
Henry Lamotte oils GmbH...	Peanut oil, refined	5592601	2
Caelo	Peanut oil, refined	5592601	3
Henry Lamotte	Peanut oil, refined	5473601	1
lamotte/Gehe	Peanut oil, refined	5592601	3
Henry Lamotte Oils GmbH	Peanut oil, refined	5697901	4
Henry Lamotte	Peanut oil, refined	5697901	2
lamotte/Phönix	Peanut oil, refined	5697901	1
Lamotte/Phoenix	Peanut oil, refined	5699601	2
Lamotte/Phoenix	Peanut oil, refined	5904100001	1
Henry Lamotte	Peanut oil, refined	5904100001	2
Henry Lamotte	Peanut oil, refined	5699601	2
Caelo	Peanut oil, refined	5699601	2
H.Lamotte/Noweda	Peanut oil, refined	5904100001	1
Henry Lamotte Oils GmbH	Peanut oil, refined	5904100091	1
Fagron	Peanut oil, refined	P101P102P103	1
Bombastus	Peanut oil, refined	264328	1
Bombastus	Peanut oil, refined	263129	1
Bombastus	Peanut oil, refined	266753	1
Bombastus	Peanut oil, refined	267966	4
Bombastus	Peanut oil, refined	271512	3
Caelo	Peanut oil, refined	7031307	1
Caelo	Peanut oil, refined	3061303	1
Henry Lamotte	Peanut oil, refined	611Erdnuss	1
Caelo	Peanut oil, refined	6061202	1
Caelo	Peanut oil, refined	466158	1
Bombastus	Peanut oil, refined	8273799	1
Caelo	Peanut oil, refined	9091303	1
Caelo	Peanut oil, refined	1013A-01954	2
Caelo	Peanut oil, refined	9101201	1
Phönix Gotha	Peanut oil, refined	101453	1
Caelo	Peanut oil, refined	10208605	1
H.Lamotte/Noweda	Peanut oil, refined	1034	1
Caelo	Peanut oil, refined	11100201	2
Caelo	Peanut oil, refined	10358602	1
Caelo	Peanut oil, refined	11163201	1
Caelo	Peanut oil, refined	11254503	1
Phönix Gotha	Peanut oil, refined	11163203	1
Caelo	Peanut oil, refined	11254504	1
Caelo	Peanut oil, refined	11313905	4
Purren Apotheke	Peanut oil, refined	11313905	1
Caelo	Peanut oil, refined	11313906	3
EuRho	Peanut oil, refined	1204018-01	1
Caelo	Peanut oil, refined	12041092	1
Noweda	Peanut oil, refined	12041902	4
Anzag am 4.05.2012	Peanut oil, refined	12041902	2
Caelo	Peanut oil, refined	12041902	12
Gehe	Peanut oil, refined	12041902	1
Caelo	Peanut oil, refined	12101502	3
Caelo	Peanut oil, refined	12080802	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Peanut oil, refined	1207211954	1
Caelo	Peanut oil, refined	12101503	10
Noweda	Peanut oil, refined	12101503	1
Erika apotheke	Peanut oil, refined	12203803	1
Caelo	Peanut oil, refined	1211503	1
Noweda	Peanut oil, refined	12203803	1
Anzag 07.12.12	Peanut oil, refined	12203803	1
Caelo	Peanut oil, refined	12274032	1
Sanacorp	Peanut oil, refined	12274302	2
Noweda	Peanut oil, refined	12274302	5
Großhandel	Peanut oil, refined	12274302	1
Phönix	Peanut oil, refined	12274302	2
Caesar/Lorenz	Peanut oil, refined	12274302	1
Caelo	Peanut oil, refined	12323204	15
AHD, 06.05.2013	Peanut oil, refined	12323204	1
Anzag	Peanut oil, refined	12323204	1
Fagron	Peanut oil, refined	12A05-N05	1
Noweda	Peanut oil, refined	123508	1
Fagron	Peanut oil, refined	12C30-N02	2
Anzag	Peanut oil, refined	12c30n02	1
Fagron	Peanut oil, refined	12c30-n02	1
Phönix Gotha	Peanut oil, refined	12C30-N02	1
Noweda	Peanut oil, refined	12C30-N02	1
Fagron	Peanut oil, refined	12H31-N05	4
Fagron	Peanut oil, refined	12J23-N01	3
Anzag	Peanut oil, refined	12H31N05	1
Fagron	Peanut oil, refined	13G30-N05	1
Fagron	Peanut oil, refined	12J23-N01	1
Caelo	Peanut oil, refined	13129202	1
Phönix Gotha	Peanut oil, refined	13129202	1
Caelo	Peanut oil, refined	13129203	41
Alliance Healthcare	Peanut oil, refined	13129203	1
Großhandel	Peanut oil, refined	13129203	1
Caelo	Peanut oil, refined	13306101	2
Gehe	Peanut oil, refined	13306103	1
Phönix Gotha	Peanut oil, refined	13306103	1
Caelo	Peanut oil, refined	13306106	19
Gehe	Peanut oil, refined	13306106	1
Phoenix	Peanut oil, refined	13306106	1
Phoenix	Peanut oil, refined	13358803	1
Anzag	Peanut oil, refined	13358804	1
Caelo	Peanut oil, refined	1335884	1
Caelo	Peanut oil, refined	1358804	1
Fagron	Peanut oil, refined	13g30-n05	1
Caesar&lorenz,Phoenix	Peanut oil, refined	14033601	1
Fagron	Peanut oil, refined	13J25B06291837	1
Caesar&lorenz,Phoenix	Peanut oil, refined	14246903	1
Caelo	Peanut oil, refined	15012016	1
Caelo	Peanut oil, refined	21121205	1
Bombastus	Peanut oil, refined	276103	2
Caelo	Peanut oil, refined	26061203	1
Caelo	Peanut oil, refined	24081208	1
Bombastus	Peanut oil, refined	277813	1
Caelo	Peanut oil, refined	3813E1954	1
Caelo	Peanut oil, refined	3112M-01954	1
Caelo	Peanut oil, refined	30101312	1
Lamotte/GEHE	Peanut oil, refined	5204402	2
Anzag 16.01.13 ek.5.13	Peanut oil, refined	5204402	2
Noweda	Peanut oil, refined	5103001	1
Fagron	Peanut oil, refined	5012A-01954	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Peanut oil, refined	5204402	1
Caelo	Peanut oil, refined	5204402	1
Henry Lamotte oils GmbH	Peanut oil, refined	5204402	2
Henry Lamotte/Jenne	Peanut oil, refined	5204402	1
Lamotte	Peanut oil, refined	5204402	2
Henry Lamotte	Peanut oil, refined	520442	1
H.Lamotte/Noweda	Peanut oil, refined	5231601	1
Fiebig 02.10.12 ek: 5,13EUR	Peanut oil, refined	5231601	1
Phönix	Peanut oil, refined	5404101	1
HenryLamotte/Phoenix	Peanut oil, refined	5231601	1
lamotte/Jenne	Peanut oil, refined	5404101	1
Lamotte / Gehe	Peanut oil, refined	5404101	1
Phönix	Peanut oil, refined	5473601	1
Lamotte	Peanut oil, refined	5473601	3
H.Lamotte/Noweda	Peanut oil, refined	5473601	1
Lamotte	Peanut oil, refined	5592601	1
Henry Lamotte	Peanut oil, refined	5592601	1
Lamotte	Peanut oil, refined	5697901	2
Lamotte	Peanut oil, refined	5699601	3
Lamotte	Peanut oil, refined	5699501	1
H.Lamotte/Noweda	Peanut oil, refined	5699601	1
Bombastus	Peanut oil, refined	80370408	1
Caelo	Peanut oil, refined	5904100001	1
Caelo	Peanut oil, refined	ch12274302	1
Caelo	Safflower oil, refined type I	16184901	1
Caelo	Safflower oil, refined type I	299143	1
Euro OTC	Soy lecithin, liquid	1604034-01	6
Euro OTC	Soy lecithin, liquid	1703001-01	8
Euro OTC	Soy lecithin, liquid	1503038-02	1
Caelo	Codliver oil (type A)	14305701	2
Caelo	Codliver oil (type A)	13396403	2
Caelo	Codliver oil (type A)	15047201	6
Lamotte/Stumpf	Codliver oil (type A)	5529601	1
Henry Lamotte Oils	Codliver oil (type A)	5529601	2
Henry Lamotte / Anzag	Codliver oil (type A)	5641901	1
Henry Lamotte Oils	Codliver oil (type A)	5693101	1
Caelo	Codliver oil (type A)	5693101	2
anzag, 9,45EUR, 19.11.14	Codliver oil (type A)	5693101	1
Caelo	Codliver oil (type A)	12346810	8
Caelo	Codliver oil (type A)	5437202	13
Caelo	Codliver oil (type A)	13106201-3392131070	1
henry lamotte oils	Codliver oil (type A)	5437202	1
Taoasis	Jojoba oil, native	137640-123934	1
Taoasis	Jojoba oil, native	133535-123399	1
Taoasis	Jojoba oil, native	141161-1819	1
Caelo	Jojoba oil, native	15116904	2
Caelo	Jojoba oil, native	15116903	5
Taoasis	Jojoba oil, native	151214-125861	2
Nachprüfung	Jojoba oil, native	14284207	2
Caelo	Jojoba oil, native	15409902	9
Caelo	Jojoba oil, native	15409903	3
Hedinger	Jojoba oil, native	15409903	1
Caelo	Jojoba oil, native	15409904	1
Caelo	Jojoba oil, native	16070601	1
Caelo	Jojoba oil, native	16070602	7
Caelo	Jojoba oil, native	16070606	5
Caelo	Jojoba oil, native	16070608	6
Caelo	Jojoba oil, native	17073309	2
Caelo	Jojoba oil, native	16070610	1
Caelo	Jojoba oil, native		2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Taoasis/Anzag	Jojoba oil, native	2a5890402-121152	1
Bombastus	Jojoba oil, native	294112	3
Bombastus	Jojoba oil, native	309480	1
Caelo	Jojoba oil, native	5756401	3
Bombastus	Jojoba oil, native	297760	2
Bombastus	Jojoba oil, native	306279	1
Lamotte	Jojoba oil, native	5756401	1
Henry Lamotte Oils GmbH	Jojoba oil, native	5913600001	2
Lamotte Sanacorp	Jojoba oil, native	8000469001	6
Lamotte Sanacorp	Jojoba oil, native	5913600001	1
Primavera	Jojoba oil, native	241H23	1
Weltecke /Noweda	Jojoba oil, native	89010416	1
Primavera Wareneingang...	Jojoba oil, native	747A23	1
Taoasis	Jojoba oil, native	133535-123190	2
Caelo	Jojoba oil, native	14052201	2
Taoasis/Anzag	Jojoba oil, native	133535-123190	1
Caelo	Jojoba oil, native	14284202	4
Fagron	Jojoba oil, native	14284202	1
Caelo	Jojoba oil, native	14284203	2
Caelo	Jojoba oil, native	14284204	6
Caelo	Jojoba oil, native	14284207	4
Caelo	Jojoba oil, native	14284208	4
Caelo	Jojoba oil, native	21041504	1
La Cura Wilco	Jojoba oil, native	2481	1
Taoasis/Anzag	Jojoba oil, native	2a5890402	1
Taoasis	Jojoba oil, native	2A5890402-121152	2
Henry Lamotte	Jojoba oil, native	5404201	1
Spinnrad	Jojoba oil, native	315601031	1
Lamotte	Jojoba oil, native	5732801	1
Lamotte/Gehe	Jojoba oil, native	5732801	1
Henry Lamotte/Alliance ...	Jojoba oil, native	5756401	1
Henry Lamotte / AHD	Jojoba oil, native	5756401	1
Taoasis/Anzag	Jojoba oil, native	5815901-118834	3
Taoasis/Anzag	Jojoba oil, native	5864501-119639	1
Henry Lamotte/ Gehe	Jojoba oil, native	5913600001	1
Caelo	Jojoba oil, native	3071301	1
Caelo	Jojoba oil, native	10071307	1
Caelo	Jojoba oil, native	11078605	2
Caelo	Jojoba oil, native	11358102	1
DAC 2006	Jojoba oil, native	11358102	1
Caelo	Jojoba oil, native	11358103	1
Caelo	Jojoba oil, native	12114701	1
Caelo	Jojoba oil, native	12114702	1
Caelo	Jojoba oil, native	12114703	1
Caelo	Jojoba oil, native	12114708	1
Caelo	Jojoba oil, native	12217003	1
Großhandel	Jojoba oil, native	12217004	1
Caelo	Jojoba oil, native	12217008	2
Caelo	Jojoba oil, native	12217009	3
Caelo	Jojoba oil, native	12217010	5
Caelo	Jojoba oil, native	13011202	5
Gehe Laborchemie Apolda...	Jojoba oil, native	13011202	1
Caelo	Jojoba oil, native	13011203	4
Caelo	Jojoba oil, native	13011204	4
Caelo	Jojoba oil, native	13129001	3
Caelo	Jojoba oil, native	13129004	5
Caelo	Jojoba oil, native	13129005	2
Caelo	Jojoba oil, native	13129006	5
Großhandel	Jojoba oil, native	13129006	1
Caelo	Jojoba oil, native	13129008	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Jojoba oil, native	15011412	1
Caelo	Jojoba oil, native	13287701	2
Caelo	Jojoba oil, native	17121302	1
Primavera	Jojoba oil, native	17009133	1
Caelo	Jojoba oil, native	19061307	1
Primavera	Jojoba oil, native	19307131	1
wilco	Jojoba oil, native	2480	1
Caelo	Jojoba oil, native	25101303	1
Bombastus	Jojoba oil, native	276323	1
Caelo	Jojoba oil, native	2613M-01986	1
HenryLamotte	Jojoba oil, native	4957101	1
HenryLamotte	Jojoba oil, native	5172601	9
Taoasis	Jojoba oil, native	4674401-67053	1
Anzag/	Jojoba oil, native	5115601	1
Lamotte/Gehe	Jojoba oil, native	5172601	1
Anzag	Jojoba oil, native	5172601	1
Henry Lamotte	Jojoba oil, native	5194301	1
HenryLamotte	Jojoba oil, native	5194301	1
Lamotte/Sanacorp	Jojoba oil, native	5194301	1
Taoasis	Jojoba oil, native	5225101-108329	1
Taoasis	Jojoba oil, native	5336301-109515	1
Taoasis	Jojoba oil, native	5336301-111368	1
Henry Lamotte	Jojoba oil, native	5429201	3
Taoasis	Jojoba oil, native	5406801-111877	1
Taoasis	Jojoba oil, native	5406801-111216	1
Taoasis GmbH	Jojoba oil, native	5445301-112580	1
Taoasis	Jojoba oil, native	5445301-112580	2
Bombastus	Jojoba oil, native	281299	1
Caeöo	Jojoba oil, native	13287702	1
Gehe/Henry Lamotte Oils...	Codliver oil (type A)	119929	1
H. Lamotte/Noweda	Codliver oil (type A)	119929	2
Henry Lamotte	Codliver oil (type A)	119929	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	123420L900	1
Henry Lamotte Oils GmbH...	Codliver oil (type A)	123823	1
Lamotte/V.d.L.	Codliver oil (type A)	123823	1
Lamotte	Codliver oil (type A)	123823	1
Henry Lamotte Oils	Codliver oil (type A)	2101204	1
Caelo	Codliver oil (type A)	3021402	1
Caelo	Codliver oil (type A)	3121201	1
Caelo	Codliver oil (type A)	5091309	1
Caelo	Codliver oil (type A)	5061202	1
Caelo	Codliver oil (type A)	6021302	1
Caelo	Codliver oil (type A)	10061302	1
Caelo	Codliver oil (type A)	11315302	2
Caelo	Codliver oil (type A)	10101207	1
Caelo	Codliver oil (type A)	11081202	1
Caelo	Codliver oil (type A)	11315303	1
Caelo	Codliver oil (type A)	12009301	1
Caelo	Codliver oil (type A)	12009302	5
Caelo	Codliver oil (type A)	12009308	2
Caelo	Codliver oil (type A)	12062501	10
Anzag 17.12.2012	Codliver oil (type A)	12062501	1
Großhandel	Codliver oil (type A)	12062501	1
Caelo	Codliver oil (type A)	12126102	1
Caelo	Codliver oil (type A)	12148703	12
Caelo	Codliver oil (type A)	121487703	1
Caelo	Codliver oil (type A)	13011301	1
Caelo	Codliver oil (type A)	13065503	4
Großhandel	Codliver oil (type A)	13106201	1
Caelo	Codliver oil (type A)	13069201	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alliance Health 21.10.2...	Codliver oil (type A)	13106201	2
Caelo	Codliver oil (type A)	13277706	8
Großhandel	Codliver oil (type A)	13277706	1
Caelo	Codliver oil (type A)	14061312	1
Caelo	Codliver oil (type A)	14075015	3
Caelo	Codliver oil (type A)	15201201	2
Caelo	Codliver oil (type A)	18031308	1
Henry Lamotte	Codliver oil (type A)	18101303	1
Henry Lamotte	Codliver oil (type A)	20021413	1
Caelo	Codliver oil (type A)	25091210	1
Caelo	Codliver oil (type A)	26061202	1
Caelo	Codliver oil (type A)	21111204	1
Caelo	Codliver oil (type A)	30041302	1
Henry Lamotte	Codliver oil (type A)	26111303	1
Henry Lamotte	Codliver oil (type A)	30041303	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	5029207	1
lamotte	Codliver oil (type A)	5029205	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	5228001	3
Henry Lamotte	Codliver oil (type A)	5029207	1
H.Lamotte/Noweda	Codliver oil (type A)	5228001	1
Henry Lamotte	Codliver oil (type A)	5228001	1
Henry Lamotte Oils	Codliver oil (type A)	5228001	1
Lamotte Jenne	Codliver oil (type A)	5228004	2
H.Lamotte/Noweda	Codliver oil (type A)	5228004	1
Henry Lamotte	Codliver oil (type A)	5228004	1
Henry Lamotte GmbH	Codliver oil (type A)	5228004	1
Hustadt Apotheke	Codliver oil (type A)	5437202	1
Henry Lamotte GmbH	Codliver oil (type A)	5437202	1
Lamotte/ Noweda	Codliver oil (type A)	5437202	7
Dr. Paulsen	Codliver oil (type A)	5437202	1
BREMER SCHLÜSSEL/ Noweda	Codliver oil (type A)	5437202	1
Henry Lamotte	Codliver oil (type A)	5437202	1
Gehe	Codliver oil (type A)	5437202	1
Henry Lamotte/V.d.L.	Codliver oil (type A)	5437202	1
Lamotte	Codliver oil (type A)	5529601	1
Henry Lamotte	Codliver oil (type A)	5529601	7
Caelo	Codliver oil (type A)	5529601	10
Bremen / Noweda	Codliver oil (type A)	5529601	1
Henry Lamotte/Phoenix	Codliver oil (type A)	5529601	2
Lamotte/V.d.L.	Codliver oil (type A)	5529601	4
Lamotte/Noweda	Codliver oil (type A)	5529601	1
henry lamotte	Codliver oil (type A)	5529601	1
Lamotte/ Noweda	Codliver oil (type A)	A458	2
Caelo	Codliver oil (type A)	2121404	1
Henry Lamotte Oils GmBh	Codliver oil (type A)	30614A	1
Caelo	Codliver oil (type A)	4061403	1
Caelo	Codliver oil (type A)	5081402	1
Henry Lamotte Oils GmBh	Codliver oil (type A)		1
Caelo	Codliver oil (type A)	8042014/A	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	8062015A	1
Caelo	Codliver oil (type A)	8071504	1
Caelo	Codliver oil (type A)	10031501	1
Caelo	Codliver oil (type A)	10061501	1
Caelo	Codliver oil (type A)	13106201	8
Caelo	Codliver oil (type A)	13277707	1
Caelo	Codliver oil (type A)	13396407	4
Caelo	Codliver oil (type A)	14305003	1
Caelo	Codliver oil (type A)	15056802	1
phönix	Codliver oil (type A)	15201201	1
Euro OTC	Codliver oil (type A)	1506009-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Codliver oil (type A)	15201203	3
Caelo	Codliver oil (type A)	15325201	2
Caesar & Loretz/Noweda	Codliver oil (type A)	15325201	1
Caelo	Codliver oil (type A)	16018202	3
phönix	Codliver oil (type A)	16074301	1
Caelo	Codliver oil (type A)	17055303	1
Caelo	Codliver oil (type A)	17011503	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	18092014A	1
Großhandel	Codliver oil (type A)	196403	1
Caelo	Codliver oil (type A)	21101401	1
Caelo	Codliver oil (type A)	23121414	1
Henry Lamotte	Codliver oil (type A)	26111402	1
lamotte/Gehe	Codliver oil (type A)	5529601	4
Caelo	Codliver oil (type A)	28041507	1
Caelo	Codliver oil (type A)	28071504	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	301214a	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	29012015A	1
pHOENIX 18.02.2014	Codliver oil (type A)	5529601	1
Henry Lamotte/Noweda	Codliver oil (type A)	5529601	1
Henry Lamotte oils GmBH...	Codliver oil (type A)	5529601	3
Noweda	Codliver oil (type A)	5529601	4
Noweda/Henry Lamotte	Codliver oil (type A)	5529601	2
lamotte/Gehe	Codliver oil (type A)	5641901	3
Henry Lamotte Oils GmBH...	Codliver oil (type A)	5641901	1
Lamotte über Sanacorp	Codliver oil (type A)	5641901	1
Henry Lamotte/Noweda	Codliver oil (type A)	5641901	1
Henry Lamotte oils GmBH...	Codliver oil (type A)	5641901	1
Henry Lamotte Oils	Codliver oil (type A)	5641901	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	5641901	1
Lamotte	Codliver oil (type A)	5693101	5
Lamotte / Gehe	Codliver oil (type A)	5693101	1
Henry Lamotte Oils/GEHE	Codliver oil (type A)	5693101	1
Lamotte/gehe	Codliver oil (type A)	5693101	3
Stumpf/Henry Lamotte	Codliver oil (type A)	5693101	1
henry lamotte	Codliver oil (type A)	5693101	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	5693101	2
Lamotte	Codliver oil (type A)	5878701	6
H.Lamotte/Noweda	Codliver oil (type A)	5878701	1
Henry Lamotte Oils/Noweda	Codliver oil (type A)	5878701	1
Stumpf / Hedinger	Codliver oil (type A)	5878701	1
Lam	Codliver oil (type A)	5878701	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	5878701	1
Caelo	Codliver oil (type A)	80000601001	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	8000601001	3
Henry Lamotte	Codliver oil (type A)	8000601001	4
Lamotte	Codliver oil (type A)	8000601001	2
lamotte jenne	Codliver oil (type A)	8000601001	1
Noweda/Henry Lamotte	Codliver oil (type A)	8000601001	2
Lamotte/Phönix	Codliver oil (type A)	8000601001	1
Caelo	Codliver oil (type A)	8000918001	6
Henry Lamotte Oils GmBH	Codliver oil (type A)	8000918001	2
Lamotte	Codliver oil (type A)	8000918001	4
Lamotte/Geibo	Codliver oil (type A)	8000918001	2
Lamotte/Anzag	Codliver oil (type A)	8000918001	1
Lamotte/Phoenix	Codliver oil (type A)	8000918001	1
Henry Lamotte	Codliver oil (type A)	8000918001	1
Henry Lamotte Oils GmBH	Codliver oil (type A)	8001176004	3
Henry Lamotte/Phoenix	Codliver oil (type A)	8001176004	1
Henry Lamotte	Codliver oil (type A)	8001176004	2
Lamotte	Codliver oil (type A)	8001176004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Codliver oil (type A)	8001176004	5
Lamotte/Gehe	Codliver oil (type A)	8001176004	1
Caelo	Codliver oil (type A)	8001215002	6
Henry Lamotte	Codliver oil (type A)	8001215002	4
Henry Lamotte Oils GmBh	Codliver oil (type A)	8001215002	1
Lamotte	Codliver oil (type A)	8001215002	1
Henry Lamotte/Noweda	Codliver oil (type A)	8001215002	1
HENRY LAMOTTE OIS GmbH	Codliver oil (type A)	800125002	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	8001765001	2
Henry Lamotte	Codliver oil (type A)	8001765001	2
Lam / Now	Codliver oil (type A)	8001765001	1
AHD	Codliver oil (type A)	8001765001	1
Caelo	Codliver oil (type A)	8001765001	2
AHD	Codliver oil (type A)	8002338002	2
Lamotte/Phoenix	Codliver oil (type A)	8002338002	5
Caelo	Codliver oil (type A)	8002338002	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	8002338002	1
Henry Lamotte/Noweda	Codliver oil (type A)	8002338002	1
Lamotte	Codliver oil (type A)	8002338002	1
Henry Lamotte Anzag	Codliver oil (type A)	8002468002	1
Lamotte/Phönix	Codliver oil (type A)	8002338002	1
Lamotte/Phönix	Codliver oil (type A)	8003727001	2
Caelo	Codliver oil (type A)	8003727001	1
Henry Lamotte	Codliver oil (type A)	800765001	1
LaMotte/Spangro	Codliver oil (type A)	80117604	1
ACM	Codliver oil (type A)	KNC0302J15	1
ACM	Codliver oil (type A)	KNC0303J15	1
ACM	Codliver oil (type A)	KNC0300J15	1
Lamotte	Codliver oil (type A)	119929	2
Phönix/Lamotte	Codliver oil (type A)		1
Stumpf/Henry Lamotte	Codliver oil (type A)	119929	1
Synopharm am 10.09.2012	Codliver oil (type A)	119929	1
Caelo	Codliver oil (type A)	119929	2
Henry Lamotte Oils GmbH	Codliver oil (type A)	123108	1
Henry Lamotte Oils, Anzag	Codliver oil (type A)	123823	2
Henry Lamotte oils GmbH	Codliver oil (type A)	123823	3
Lamotte/Phönix	Codliver oil (type A)	123823	1
Henry Lamotte/Sanacorp	Codliver oil (type A)	123823	1
Lamotte/L3	Codliver oil (type A)	123823	1
Stumpf/Henry Lamotte	Codliver oil (type A)	123823	1
lamotte/Gehe	Codliver oil (type A)	123823	1
Phoenix	Codliver oil (type A)	11315303	1
Caelo	Codliver oil (type A)	8131-01985	1
Henry Lamotte Oils GmBh	Codliver oil (type A)	6082013A	1
Fagron	Codliver oil (type A)	12H15-N02	1
Caelo	Codliver oil (type A)	13347501	1
Henry Lamotte Oils GmBh	Codliver oil (type A)	16012014A	1
Gehe	Codliver oil (type A)	2613A-01985	1
Henry Lamotte Oils GmBh	Codliver oil (type A)	280913A	1
Henry Lamotte Oils GmBh	Codliver oil (type A)	291013A	1
Noweda	Codliver oil (type A)	4537202	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	4931805	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	5029205	3
Lamotte/ Sanacorp	Codliver oil (type A)	5029207	1
Noweda/Lamotte Oils	Codliver oil (type A)	5029207	1
Caelo	Codliver oil (type A)	5228001	2
Noweda	Codliver oil (type A)	5228001	1
Lamotte Oils GmbH / G	Codliver oil (type A)	5228004	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	5228004	2
Lamotte	Codliver oil (type A)	5228004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Codliver oil (type A)	5228004	2
Lamotte/GEHE	Codliver oil (type A)	5437202	8
Henry Lamotte Oils GmbH	Codliver oil (type A)	5437202	3
Noweda / Henry Lamotte	Codliver oil (type A)	5437202	1
Lamotte/gehe	Codliver oil (type A)	5437202	4
Lamotte Oils GmbH / E	Codliver oil (type A)	5437202	2
Lamotte	Codliver oil (type A)	5437202	4
Henry Lamotte/Sanacorp	Codliver oil (type A)	5529601	1
Henry Lamotte/Kehr	Codliver oil (type A)	5529601	1
Henry Lamotte Oils GmbH	Codliver oil (type A)	5529601	1
Henry Lamotte Oils/Noweda	Codliver oil (type A)	5529601	1
Hedinger	Codliver oil (type A)	8000601001	1
Henry Lamotte Oils, Gehe	Codliver oil (type A)	5529601	1
henry lamotte oils gmbh	Codliver oil (type A)	8000918001	1
henry lamotte oils gmbh	Codliver oil (type A)	8001176004	1
EUro/AHD	Soy lecithin, liquid	1305019-01	1
Gehe	Soy lecithin, liquid	21022111/2	1
Lamotte/Noweda	Flaxseed oil, native	103229	1
Henry/Jenne	Flaxseed oil, native	103277	3
Lamotte	Flaxseed oil, native	124116	2
Euro OTC	Flaxseed oil, native	124116	1
Henry Lamotte/Noweda	Flaxseed oil, native	124116	1
Henry/AHD	Flaxseed oil, native	124337	2
Caelo	Flaxseed oil, native	9062015A	1
Caelo	Flaxseed oil, native	13361202	4
Caelo	Flaxseed oil, native	14087903	1
Bombastus	Flaxseed oil, native	14087903	1
Caelo	Flaxseed oil, native	14087908	4
Caelo	Flaxseed oil, native	14087909	1
Caelo	Flaxseed oil, native	14306808	1
Caelo	Flaxseed oil, native	14306811	2
Caelo	Flaxseed oil, native	15116702	5
Caelo	Flaxseed oil, native	15175608	4
Bombastus	Flaxseed oil, native	285121	1
Gehe	Flaxseed oil, native	288025	1
Bombastus	Flaxseed oil, native	288025	1
Caesar & Loretz GmbH Hi...	Flaxseed oil, native	289164	1
Bombastus	Flaxseed oil, native	291126	1
Bombastus	Flaxseed oil, native	291900	1
Bombastus	Flaxseed oil, native	293950	1
statfold	Flaxseed oil, native	29409	1
Bombastus	Flaxseed oil, native	297286	1
statfold	Flaxseed oil, native	33588	2
Lamotte/Otto	Flaxseed oil, native	5757301	1
ACM	Flaxseed oil, native	KNC0446J15	1
Henry Lamotte	Flaxseed oil, native	1000003255	1
Henry Lamotte	Flaxseed oil, native	1000003598	1
Henry Lamotte Oils GmbH	Flaxseed oil, native	1000005198	1
Caelo	Flaxseed oil, native	15415906	5
Gehe/Vitalhaus	Flaxseed oil, native	160702	1
Caelo	Flaxseed oil, native	16321106	4
Caelo	Flaxseed oil, native	17046711	1
Bombastus	Flaxseed oil, native	17002089	1
Bombastus	Flaxseed oil, native	296730	1
Bombastus	Flaxseed oil, native	302269	2
Bombastus	Flaxseed oil, native	303096	3
Bombastus	Flaxseed oil, native	303440	1
Bombastus	Flaxseed oil, native	304454	2
Bombastus	Flaxseed oil, native	306179	2
Bombastus	Flaxseed oil, native	266111	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Flaxseed oil, native	307988	1
Bombastus	Flaxseed oil, native	308634	1
Bombastus	Flaxseed oil, native	271645	2
Bombastus	Flaxseed oil, native	272528	1
Kruger/Lamotte/30.12.13	Flaxseed oil, native	102503	1
Noweda / Henry Lamotte	Flaxseed oil, native	119529	1
Gehe	Flaxseed oil, native	273676	1
Noweda / Henry Lamotte	Flaxseed oil, native	119677	2
Phönix 11.09.2012	Flaxseed oil, native	118553	1
Stumpf/Henry Lamotte	Flaxseed oil, native	124116	1
Caelo	Flaxseed oil, native	12217504	2
Caelo	Flaxseed oil, native		2
Caelo	Flaxseed oil, native	7111204	1
Gehe	Flaxseed oil, native	12041701	1
Caelo	Flaxseed oil, native	12041701	1
Caelo	Flaxseed oil, native	12217505	2
Caelo	Flaxseed oil, native	12274205	2
Henry Lamotte Oils GMBH...	Flaxseed oil, native	12D05-N02	1
Gehe	Flaxseed oil, native	13178203	1
Caelo	Flaxseed oil, native	13178203	4
Caelo	Flaxseed oil, native	13178209	3
Caelo	Flaxseed oil, native	13361203	2
Caelo	Flaxseed oil, native	23021302	1
Bombastus	Flaxseed oil, native	275487	2
Bombastus	Flaxseed oil, native	276106	1
Bombastus	Flaxseed oil, native	276727	2
Bombastus	Flaxseed oil, native	276105	1
Bombastus	Flaxseed oil, native	279769	1
Bombastus	Flaxseed oil, native	280711	1
Bombastus	Flaxseed oil, native	281489	2
s.Prüfprotokoll	Flaxseed oil, native		1
Bomb/ Noweda	Flaxseed oil, native	290569	1
Bombastus	Flaxseed oil, native	282391	1
Bombastus	Flaxseed oil, native	274875	1
Caelo	Laurel oil	14176301	2
Caelo	Laurel oil	14176303	1
Caelo	Laurel oil	16355701	1
Caesar & Loretz/Noweda	Laurel oil	15216401	1
Caelo	Laurel oil	16355702	1
Caelo	Corn oil, refined	15296407	4
Caelo	Corn oil, refined	15050506	1
Caelo	Corn oil, refined	11267907	1
Resana	Corn oil, refined	1208-3678	1
Caelo	Corn oil, refined	12141511	2
Caelo	Corn oil, refined	13010707	2
Caelo	Corn oil, refined	13010710	1
Resana	Corn oil, refined	1308-4063	1
Primavera WE 29.03.17 -...	Almond oil, refined	241C25	1
Henry Lamotte/Anzag	Almond oil, refined	416E-01950	1
Bahnhof-Apotheke	Almond oil, refined	-30-160809	1
Henry Lamotte Oils / Gehe	Almond oil, refined	1000002454	1
Henry Lamotte	Almond oil, refined	1000002454	1
Henry Lamotte/ Henry La...	Almond oil, refined	1000003242	19
Henry Lamotte/Sanacorp	Almond oil, refined	1000002454	2
Lamotte	Almond oil, refined	1000003242	2
Henry Lamotte/Gehe	Almond oil, refined	1000003242	1
Euro OTC	Almond oil, refined	1000003242	2
WE:Phönix 25.08.2016 EK...	Almond oil, refined	1000003242	2
Henry Lamotte Oils GmbH	Almond oil, refined	1000003242	2
Lamotte / Noweda	Almond oil, refined	1000003242	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils / Gehe	Almond oil, refined	1000003242	1
Henry Lamotte	Almond oil, refined	1000003242	3
Henry Lamotte/Noweda	Almond oil, refined	1000003242	1
Lamotte/Sanacorp	Almond oil, refined	1000003242	2
L20/Lamotte	Almond oil, refined	1000003242	1
Henry Lamotte Oils GmbH...	Almond oil, refined	1000004431	2
L 20/ Henry Lamotte	Almond oil, refined	1000004431	1
Gehe/Henry Lamotte	Almond oil, refined	1000004431	1
Caelo	Almond oil, refined	1000004431	2
Henry Lamotte Oils GmbH	Almond oil, refined	1000004431	1
Henry Lamotte/ Noweda	Almond oil, refined	1000004431	1
Henry Lamotte Oils	Almond oil, refined	1000004431	2
Hedinger	Almond oil, refined	1000004431	1
Henry Lamotte/Noewda	Almond oil, refined	1000004431	1
Lamotte	Almond oil, refined	1000004815	1
Spinnrad	Almond oil, refined	1215048093	1
Euro OTC	Almond oil, refined	1510014-01	3
Fagron	Almond oil, refined	12345	1
Caelo	Almond oil, refined	15197401	10
Caelo	Almond oil, refined	15197403	10
Fagron	Almond oil, refined	15197404	1
Caelo	Almond oil, refined	15197406	6
Caelo	Almond oil, refined	15197404	4
Caelo	Almond oil, refined	15197407	8
Caelo	Almond oil, refined	15197409	2
Bombastus	Almond oil, refined	15197412	1
Caelo	Almond oil, refined	15331611	1
Caelo	Almond oil, refined	15409105	1
Caelo	Almond oil, refined	15409701	4
Fagron	Almond oil, refined	15409702	1
Caelo	Almond oil, refined	15409702	4
Caelo	Almond oil, refined	15409703	7
Caeolo/Noweda	Almond oil, refined	15409703	1
Caelo	Almond oil, refined	15409704	9
Alliance H.	Almond oil, refined	15409704	1
Caelo	Almond oil, refined	15409705	9
Fagron	Almond oil, refined	15409705	1
Caelo	Almond oil, refined	15409707	1
Caelo	Almond oil, refined	15409710	10
Phönix	Almond oil, refined	15409710	1
Caelo	Almond oil, refined	15409713	5
Caelo	Almond oil, refined	15409714	6
Fagron	Almond oil, refined	15H03-B07-314457	3
Euro OTC	Almond oil, refined	1603025-01	3
Caelo	Almond oil, refined	16071001	2
Euro OTC	Almond oil, refined	16071002	1
Caelo	Almond oil, refined	16071002	11
Caelo	Almond oil, refined	16071003	11
Caelo	Almond oil, refined	16071004	5
Caelo	Almond oil, refined	16071005	3
Caelo	Almond oil, refined	16071006	5
Caelo	Almond oil, refined	16071007	6
Cealo	Almond oil, refined	16071007	3
Alte Apotheke	Almond oil, refined	16071007	1
Caelo	Almond oil, refined	16071008	1
Caelo	Almond oil, refined	16071009	1
Caelo	Almond oil, refined	16144403	4
Caelo	Almond oil, refined	16327302	16
Caelo	Almond oil, refined	16149802	1
Caelo	Almond oil, refined	16327303	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Almond oil, refined	16327304	8
Caelo	Almond oil, refined	16327305	6
Caelo	Almond oil, refined	16327308	4
Caesar & Loretz GmbH	Almond oil, refined	16327308	1
Caelo	Almond oil, refined	16327311	4
Caelo	Almond oil, refined	1633327302	1
Caelo	Almond oil, refined	1637302	1
Fagron	Almond oil, refined	16A05-B02-320278	1
Henry Lamotte Oils GMBH	Almond oil, refined	170316M8001497001	1
Caelo	Almond oil, refined	17042203	8
Caelo	Almond oil, refined	17042204	4
Caelo	Almond oil, refined	1716506	1
Caelo	Almond oil, refined	17176506	2
Caelo	Almond oil, refined	17176501	3
Beiersdorf	Almond oil, refined	17207905	1
Caelo	Almond oil, refined	17207908	1
Caelo	Almond oil, refined	17176511	1
Caelo	Almond oil, refined	18015007	3
Caelo	Almond oil, refined	18070604	1
Bombastus	Almond oil, refined	291509	5
Noweda	Almond oil, refined	294784	1
Caelo	Almond oil, refined	295616	3
Bombastus	Almond oil, refined	297872	3
Bombastus	Almond oil, refined	298744	2
Bombastus	Almond oil, refined	298745	7
Bombastus	Almond oil, refined	300841	2
Bombastus	Almond oil, refined	300840	1
Bombastus	Almond oil, refined	302210	3
Bombastus	Almond oil, refined	302211	3
Bombastus	Almond oil, refined	306291	1
Bombastus	Almond oil, refined	304023	1
Bombastus	Almond oil, refined	307298	2
Bombastus	Almond oil, refined	306292	1
Bombastus	Almond oil, refined	307299	2
Bombastus	Almond oil, refined	308919	1
Bombastus	Almond oil, refined	310859	1
Taoasis	Almond oil, refined	31952-123336	1
Taoasis	Almond oil, refined	33489-124674	1
Lamotte/Phönix	Almond oil, refined	3243	1
Taoasis	Almond oil, refined	34371-125478	1
Lamotte/Phoenix	Almond oil, refined	8000712001	4
Caelo	Almond oil, refined	71/AKD15	1
Caelo	Almond oil, refined	72/AKD16	1
Taoasis	Almond oil, refined	35171-125735	1
Henry Lamotte	Almond oil, refined	8000712001	8
Lamotte/Gehe	Almond oil, refined	8000725001	1
Henry Lamotte Oils / Gehe	Almond oil, refined	8000725001	2
Henry Lamotte Oils GmbH	Almond oil, refined	8000725001	3
Henry Lamotte Oils GmbH	Almond oil, refined	8000725001	1
Henry Lamotte	Almond oil, refined	8000725001	2
Lamotte über Sanacorp	Almond oil, refined	8000915001	3
Caelo	Almond oil, refined	8000915001	5
Henry Lamotte/ Henry La...	Almond oil, refined	8000915001	16
Henry Lamotte Oils GmbH	Almond oil, refined	8000915001	3
Henry Lamotte	Almond oil, refined	8000915001	4
henry lamotte oils Gmh/N	Almond oil, refined	8000915001	2
Henry Lamotte / Allianc...	Almond oil, refined	8000915001	2
Noweda/Lamotte	Almond oil, refined	8000915001	1
Lamotte	Almond oil, refined	8000915001	2
Henry Lamotte/ Noweda	Almond oil, refined	8000915001	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils / Gehe	Almond oil, refined	8000915002	2
Henry Lamotte Oil	Almond oil, refined	8000915001	1
henry lamotte oils Gmh/N	Almond oil, refined	8000915002	1
Lamotte über Sanacorp	Almond oil, refined	8000915002	1
Noweda/Henry Lamotte	Almond oil, refined	8000915001	1
Henry Lamotte/Sanacorp	Almond oil, refined	8000915002	8
Caelo	Almond oil, refined	8000915002	5
Henry Lamotte	Almond oil, refined	8000915002	3
Lamotte/Gehe	Almond oil, refined	8000915002	2
Henry Lamotte Oils GmbH	Almond oil, refined	8000915002	4
Henry Lamotte Oil	Almond oil, refined	8000915002	1
Lamotte	Almond oil, refined	8000915002	2
Henry Lamotte Olis GmbH	Almond oil, refined	8000915002	2
Henry Lamotte Oils GmbH	Almond oil, refined	8000915002	1
Gehe/henry	Almond oil, refined	8000915002	1
Henry Lamotte	Almond oil, refined	8001497001	3
Lamotte/ Noweda	Almond oil, refined	8001497001	2
Anzag, Henry Lamotte	Almond oil, refined	8001497001	1
henry lamotte	Almond oil, refined	8001497001	1
Lamotte/Gehe	Almond oil, refined	8001497001	2
Henry Lamotte Oils / Gehe	Almond oil, refined	8001497001	1
Henry Lamotte Oils GmbH	Almond oil, refined	8001497001	1
Henry Lamotte/Sanacorp	Almond oil, refined	8001497001	3
Lamotte	Almond oil, refined	8001497001	1
Henry Lamotte Sanacorp	Almond oil, refined	8001497001	1
Sanacorp WE:04.02.2017 ...	Almond oil, refined	8001497001	1
Henry Lamotte Oils GmbH	Almond oil, refined	8001497001	2
Henry Lamotte / Noweda	Almond oil, refined	8001497001	2
Henry Lamotte Oils	Almond oil, refined	8001497001	1
lamotte/gehe	Almond oil, refined	8001497001	2
Henry Lamotte Oils GmbH...	Almond oil, refined	8002653001	1
Bombastus	Almond oil, refined	8002653001	2
Henry Lamotte Oils GmbH	Almond oil, refined	8002653001	1
Henry Lamotte Oils / Gehe	Almond oil, refined	8002653001	3
L 20/ Henry Lamotte	Almond oil, refined	8002653001	2
Caelo	Almond oil, refined	8002653001	3
Finteler Apotheke	Almond oil, refined	8002653001	1
lamotte/gehe	Almond oil, refined	8002653001	1
Henry Lamotte	Almond oil, refined	8002653001	2
L20/ Henry Lamotte	Almond oil, refined	8002675001	1
Henry Lamotte/Noweda	Almond oil, refined	8002653001	1
Henry Lamotte Oils	Almond oil, refined	8002675001	2
Lamotte/Noweda	Almond oil, refined	8002675001	1
Lamotte	Almond oil, refined	8002675001	4
henry lamotte oils	Almond oil, refined	8002675001	1
Henry Lamotte	Almond oil, refined	8003167001	2
Sanacorp/Henry Lamotte ...	Almond oil, refined	8003167001	3
L10/HENRY LAMOTTE	Almond oil, refined	8003167001	1
Henry Lamotte Oils GmbH...	Almond oil, refined	8003167001	1
Henry Lamotte/Noewda	Almond oil, refined	8003503001	2
Lamotte / Noweda	Almond oil, refined	8003167001	1
Lamotte	Almond oil, refined	8003503001	1
henry lamotte/gehe	Almond oil, refined	8003503001	1
Lamotte/Sanacorp	Almond oil, refined	8003503002	2
henry lamotte oils	Almond oil, refined	8003503001	1
Sanacoro W.E:16.03.18 ...	Almond oil, refined	8003503002	1
Henry Lamotte Oils / Gehe	Almond oil, refined	8003503002	2
Lamotte/Gehe	Almond oil, refined	8003503002	1
L20/ Henry Lamotte	Almond oil, refined	8003503002	1
EuRho/Noweda	Almond oil, refined	8003885002	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte/ Alliance	Almond oil, refined	8003885002	1
BioPräp	Almond oil, refined	E1610090	1
ACM	Almond oil, refined	KNC1560J16	1
Primavera	Almond oil, refined	11G23	2
Primavera	Almond oil, refined	82B23	2
Primavera	Almond oil, refined	218D24	1
Primavera	Almond oil, refined	21D24	1
Caelo	Almond oil, refined	465940	1
Primavera	Almond oil, refined	480K23	1
Primavera Wareneingang...	Almond oil, refined	792D23	1
Primavera 19.10.15	Almond oil, refined	944H23	1
Taoasis	Almond oil, refined	12842-116104	1
Eigenherstellung	Almond oil, refined	20915HML01	1
Caelo	Almond oil, refined	40315PML01	1
Kehr / Lamotte	Almond oil, refined	50115PML01	1
Eigenherstellung	Almond oil, refined	50115PML02	1
Caesar & Loretz GmbH Hi...	Almond oil, refined	510014-01	1
Caelo	Almond oil, refined	8041501	1
Primavera	Almond oil, refined	7708L22	1
Eig	Almond oil, refined	80515HML01	1
Prima Vera	Almond oil, refined	8603152	1
Euro OTC	Almond oil, refined	540014-01	1
H. Lamotte/Noweda	Almond oil, refined	1000002454	1
Lamotte/AHD	Almond oil, refined	1000003242	2
Eigenherstellung	Almond oil, refined	120914HAS01	1
Caelo	Almond oil, refined	10061502	1
Eigenherstellung	Almond oil, refined	120814HAS01	1
Caelo	Almond oil, refined	123test	1
Henry Lamotte Oils GmBh	Almond oil, refined	13032014B	1
Taoasis	Almond oil, refined	130973-116149	1
Euro OTC	Almond oil, refined	1310022-01	1
Taoasis	Almond oil, refined	130973-118607	1
Euro OTC	Almond oil, refined	131002202	1
Caelo	Almond oil, refined	13236411	7
Aliance Hearthcave 11.0...	Almond oil, refined	13236402	1
Taoasis	Almond oil, refined	132805-118635	1
Caelo	Almond oil, refined	13326411	1
Fagron	Almond oil, refined	13350001	2
Caelo	Almond oil, refined	13350002	12
Ph.Eur.7.0	Almond oil, refined	13350003	1
Caelo	Almond oil, refined	13350006	7
Taoasis	Almond oil, refined	133805-120282	1
Caelo	Almond oil, refined	13350007	4
Taoasis	Almond oil, refined	134032-120671	1
Taoasis	Almond oil, refined	134475-1215919	1
Caelo	Almond oil, refined	1350003	1
Taoasis	Almond oil, refined	135831-121621	2
Fagron	Almond oil, refined	13K28-B01-288270	1
Caelo	Almond oil, refined	14042702	3
Caelo/ Noweda	Almond oil, refined	14042704	1
Caelo	Almond oil, refined	14042705	7
Caelo	Almond oil, refined	14042704	3
Caelo	Almond oil, refined	14042706	8
Bombastus	Almond oil, refined	14042706	1
Caelo	Almond oil, refined	14042707	5
Caelo	Almond oil, refined	14042708	7
Caelo	Almond oil, refined	14042717	3
Caelo	Almond oil, refined	14042714	1
Euro OTC	Almond oil, refined	1407050-01	1
Euro OTC	Almond oil, refined	1407050-02	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Almond oil, refined	14139206	1
Caelo	Almond oil, refined	14139304	2
AHD	Almond oil, refined	14139305	1
Caelo	Almond oil, refined	14139306	2
Caelo	Almond oil, refined	14139309	3
Apolda/ Sanacorp	Almond oil, refined	14139309	1
Sanacorp	Almond oil, refined	1413935	1
Caelo	Almond oil, refined	1421508	1
Caelo	Almond oil, refined	14215804	4
Caesar & Loretz GmbH	Almond oil, refined	14215804	1
Caelo	Almond oil, refined	14215805	9
Caelo	Almond oil, refined	14215808	6
Apolda/ Sanacorp	Almond oil, refined	14215807	1
Gehe 22.12.14	Almond oil, refined	14215808	1
Caelo	Almond oil, refined	14215809	10
Fagron	Almond oil, refined	14215808	1
Gehe Kassel	Almond oil, refined	14215814	1
Caelo	Almond oil, refined	14215814	4
Caelo	Almond oil, refined	14215815	2
Caelo	Almond oil, refined	14215817	6
Fagron	Almond oil, refined	14D24-B07	1
Caelo	Almond oil, refined	1439310	1
Fagron	Almond oil, refined	14D24-B07-296709	2
Fagron	Almond oil, refined	14F04-B01	2
Fagron	Almond oil, refined	14F04-B01-294721	2
Fagron	Almond oil, refined	14J01-B03	1
Fagron	Almond oil, refined	14J01-B03-302327	2
Caelo	Almond oil, refined	15071406	1
Sanacorp	Almond oil, refined	15078401	1
Caelo	Almond oil, refined	15078402	4
Caelo	Almond oil, refined	15078403	8
Holdermann	Almond oil, refined	15078403	1
GEHE 11.05.	Almond oil, refined	15078404	1
Caelo	Almond oil, refined	15078405	9
Sanacorp	Almond oil, refined	15078405	1
Caelo	Almond oil, refined	15078406	3
Euro	Almond oil, refined	1510014-01	1
Caelo	Almond oil, refined	15078407	2
Caelo	Almond oil, refined	15127302	2
Fagron	Almond oil, refined	15127303	1
Caelo	Almond oil, refined	15127304	6
Caelo	Almond oil, refined	15127305	6
caesar&loretz,Phönix,11...	Almond oil, refined	15127304	1
Caelo	Almond oil, refined	15127306	13
Caelo	Almond oil, refined	15127307	4
Caelo	Almond oil, refined	15197402	8
Caelo	Almond oil, refined	15197405	2
Caelo	Almond oil, refined	15197408	6
Caelo	Almond oil, refined	15409706	1
Fagron	Almond oil, refined	15F04-B01	1
Fagron	Almond oil, refined	15F04-B01-309291	3
Caelo	Almond oil, refined	16021510	1
Caelo	Almond oil, refined	17061403	1
Henry Lamotte Oils GmbH	Almond oil, refined	17082015B	1
Caelo	Almond oil, refined	17111402	1
Caelo	Almond oil, refined	18081402	1
Eigenherstellung	Almond oil, refined	180815HJK01	1
Eigenherstellung	Almond oil, refined	180915HML01	1
Caelo	Almond oil, refined	220515PML01	1
Caelo	Almond oil, refined	230315HML01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Eigenherstellung	Almond oil, refined	260515HML01	1
Henry Lamotte Oils GmBH	Almond oil, refined	260814A	1
Bombastus	Almond oil, refined	283801	2
Bombastus	Almond oil, refined	284968	2
Bombastus	Almond oil, refined	286254	5
Caelo	Almond oil, refined	284968	1
Bombastus	Almond oil, refined	288428	3
Bombastus	Almond oil, refined	289709	1
Bombastus	Almond oil, refined	290136	1
Eigenherstellung	Almond oil, refined	290714HJK01	1
Eigenherstellung	Almond oil, refined	290714HJK02	1
Caelo	Almond oil, refined	29081404	1
Bombastus	Almond oil, refined	292591	2
Bombastus	Almond oil, refined	292592	6
Bombastus	Almond oil, refined	294784	3
Bombastus	Almond oil, refined	295616	2
Bombastus	Almond oil, refined	296519	2
Caelo	Almond oil, refined	31031403	1
Caelo	Almond oil, refined	30121416	1
Bombastus	Almond oil, refined	47/2014	1
Henry Lamotte Oils GmbH	Almond oil, refined	5415901	2
Bombastus	Almond oil, refined	5415901	3
Henry Lamotte / Noweda	Almond oil, refined	5415901	1
Lamotte/Gehe	Almond oil, refined	5415901	7
Caelo	Almond oil, refined	5415901	9
Großhandel	Almond oil, refined	5415901	2
Henry Lamotte Oils	Almond oil, refined	5415901	2
Henry Lamotte Oils GmbH...	Almond oil, refined	5415901	1
Lamotte	Almond oil, refined	5415901	9
Henry Lamotte	Almond oil, refined	5615001	7
Lamotte	Almond oil, refined	561001	1
Caelo	Almond oil, refined	5615001	7
Gehe/Henry Lamotte Oils...	Almond oil, refined	5615001	2
Henry Lamotte Oils GmbH...	Almond oil, refined	5615001	1
Lamotte / Alliance Heal...	Almond oil, refined	5615001	2
Lamotte / Gehe	Almond oil, refined	5615001	1
lamotte/Gehe	Almond oil, refined	5615001	1
Lamotte/Sanacorp	Almond oil, refined	5615001	2
Lamotte/ Noweda	Almond oil, refined	5615001	3
Phönix/ Lamotte	Almond oil, refined	5615001	1
Lamotte / Noweda	Almond oil, refined	5615001	1
H.Lamotte/Noweda	Almond oil, refined	5615001	3
Henry Lamotte Oils GmbH...	Almond oil, refined	5615001	1
Lamotte/Phönix	Almond oil, refined	5615003	1
LAMOTTE/NOWEDA	Almond oil, refined	5615003	4
Lamotte/Gehe	Almond oil, refined	5615003	2
Lamotte/Phoenix	Almond oil, refined	5615003	2
Phönix/ Lamotte	Almond oil, refined	5615003	2
Noweda/Henry Lamotte	Almond oil, refined	5615003	4
Phönix/Lamotte	Almond oil, refined	5615003	5
lamotte	Almond oil, refined	5615003	1
Henry Lamotte Oils GmbH...	Almond oil, refined	5615003	1
Lamotte	Almond oil, refined	5615003	4
Henry Lamotte Oils GmbH	Almond oil, refined	5615003	1
Henry Lamotte Oils / Gehe	Almond oil, refined	5615003	2
Henry Lamotte	Almond oil, refined	5615003	3
Lamotte/Noweda	Almond oil, refined	5615003	2
H. Lamotte/Noweda	Almond oil, refined	5615003	1
Henry Lamotte	Almond oil, refined	5669601	1
Caelo	Almond oil, refined	5669601	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte/ Sana	Almond oil, refined	5669601	1
Lamotte	Almond oil, refined	5669601	1
henry lamotte	Almond oil, refined	5675003	1
Henry Lamotte Oils GmbH	Almond oil, refined	5703001	2
Caelo	Almond oil, refined	5703001	2
Henry Lamotte Olis GmbH...	Almond oil, refined	5703001	6
Lamotte/Phoenix	Almond oil, refined	5703001	1
Lamotte/Sanacorp	Almond oil, refined	5703001	1
Henry Lamotte / Noweda	Almond oil, refined	5703001	1
Lamotte über Phoenix	Almond oil, refined	5703001	1
Lamotte/Sanacorp	Almond oil, refined	5703101	1
Henry Lamotte	Almond oil, refined	5703101	3
Lamotte	Almond oil, refined	5703101	2
Lamotte/Gehe	Almond oil, refined	5703101	2
Euro OTC	Almond oil, refined	5703101	1
Lamotte / Gehe	Almond oil, refined	5703101	2
Lamotte	Almond oil, refined	5713001	1
Henry Lamotte	Almond oil, refined	5748401	1
Noweda/Henry Lamotte	Almond oil, refined	5748501	1
Stumpf/Henry Lamotte	Almond oil, refined	5748601	2
Lamotte/Sanacorp	Almond oil, refined	5748601	6
Lamotte / Noweda	Almond oil, refined	5748601	2
Henry Lamotte Oils GmbH...	Almond oil, refined	5748601	3
Lamotte	Almond oil, refined	5748601	7
Noweda / Henry Lamotte	Almond oil, refined	5748601	1
Henry Lamotte	Almond oil, refined	5748601	7
H.Lamotte/Noweda	Almond oil, refined	5748601	2
Henry Lamotte Oils	Almond oil, refined	5748601	2
Lamotte / Alliance Heal...	Almond oil, refined	5748601	1
Caelo	Almond oil, refined	5748601	2
Lamotte über Phoenix	Almond oil, refined	5748601	3
Noweda/Lamotte	Almond oil, refined	5748601	1
Lamotte/Phoenix	Almond oil, refined	5748601	2
Henry La Motte	Almond oil, refined	5748601	2
lamotte/Phönix	Almond oil, refined	5748601	1
Lamotte/ Phoenix	Almond oil, refined	5748601	1
henry lamotte	Almond oil, refined	5748601	1
Henry Lamotte/Alliance H	Almond oil, refined	5748601	1
Noweda	Almond oil, refined	5748601	2
Lamotte/Gehe	Almond oil, refined	5748601	2
Henry Lamotte Oils/Sana...	Almond oil, refined	5919600001	3
Stumpf/Henry Lamotte	Almond oil, refined	5919600001	2
Henry Lamotte	Almond oil, refined	5919600001	10
Henry Lamotte Olis GmbH...	Almond oil, refined	5919600001	1
Lamotte	Almond oil, refined	5919600001	6
Henry Lamotte Oils	Almond oil, refined	5919600001	1
Henry Lamotte Oils / Gehe	Almond oil, refined	5919600001	3
Lamotte/Sanacorp 02.11...	Almond oil, refined	5919600001	1
lamotte/phönix	Almond oil, refined	5919600001	1
Noweda/Lamotte	Almond oil, refined	5919600001	1
Henry Lamotte/ Noweda	Almond oil, refined	5919600001	1
Noweda Lamotte GmbH	Almond oil, refined	5919600001	1
Henry Lamotte Oils GmbH	Almond oil, refined	5919600001	1
H.Lamotte/Noweda	Almond oil, refined	5919600001	1
Lamotte/Anzag	Almond oil, refined	8000712001	4
Lamotte/Sanacorp	Almond oil, refined	8000712001	4
Lamotte	Almond oil, refined	8000712001	10
Lamotte/Noweda	Almond oil, refined	8000712001	2
AHD / H.Lamotte	Almond oil, refined	8000712001	1
Lamotte/PPL	Almond oil, refined	8000712001	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda/Lamotte	Almond oil, refined	8000712001	3
Henry lamotte	Almond oil, refined	8000712001	1
lamotte/noweda	Almond oil, refined	8000712001	2
Lamotte über Phoenix	Almond oil, refined	8000712001	2
Henry Lamotte Oils GmbH	Almond oil, refined	8000712001	1
henry lamotte	Almond oil, refined	8000712001	3
Caelo	Almond oil, refined	8000712001	1
Henry Lamotte Oils / No...	Almond oil, refined	8000712001	1
Henry Lamotte/ Gehe	Almond oil, refined	8000712001	1
Henry Lamotte / Kehr	Almond oil, refined	8000712001	2
Pharmachem/Noweda	Almond oil, refined	8000712001	1
H.Lamotte/Noweda	Almond oil, refined	8000712001	1
Lamotte/Gehe	Almond oil, refined	8000712001	1
Henry Lamotte/Noweda	Almond oil, refined	8000712001	1
Noweda / Henry Lamotte	Almond oil, refined	8000725001	1
Henry Lamotte/Kehr	Almond oil, refined	8000725001	1
Lamotte	Almond oil, refined	8000725001	2
Henry Lamotte/Noweda	Almond oil, refined	8000725001	1
Henry Lamotte Oils	Almond oil, refined	8000725001	1
Henry Lamotte Olis GmbH...	Almond oil, refined	8000915001	2
henry lamotte	Almond oil, refined	8000915001	1
Henry Lamotte Oils	Almond oil, refined	8000915001	1
AHD/ Henry Lamotte	Almond oil, refined	8000915002	1
phönix	Almond oil, refined	8001497001	1
Lamotte	Almond oil, refined	8002653001	1
Lamotte/Sanacorp	Almond oil, refined	8748401	1
Biopräp / GEHE	Almond oil, refined	E1403027	1
Biopräp / GEHE	Almond oil, refined	E1410088	1
Caelo	Almond oil, refined		2
ACM	Almond oil, refined	KNC0766J15	1
ACM	Almond oil, refined	KNC1364J15	1
ACM	Almond oil, refined	KNC1100J15	1
Bombastus	Almond oil, refined	264164	1
Bombastus	Almond oil, refined	267099	2
Bombastus	Almond oil, refined	269484	1
Bombastus	Almond oil, refined	270829	1
Bombastus	Almond oil, refined	270830	1
Bombastus	Almond oil, refined	273578	4
anzag, 05.11.2012	Almond oil, refined	273578	1
Phönix 15.03.2013	Almond oil, refined	273579	1
Bombastus	Almond oil, refined	274564	1
Taoasis	Almond oil, refined	102235-111893	1
HenryLamotte/Gehe	Almond oil, refined	130145	3
Henry Lamotte GmbH, Mer...	Almond oil, refined	130145	1
Fiebig 21.02.13 ek: 2,16EUR	Almond oil, refined	130145	1
Henry Lamotte Oils GmbH	Almond oil, refined	130145	2
Phönix 16.03.2013	Almond oil, refined	130145	2
Caelo	Almond oil, refined	130145	4
Henry Lamotte	Almond oil, refined	130145	3
Lamotte/Noweda	Almond oil, refined	130145	3
H.Lamotte/Noweda	Almond oil, refined	130145	2
Henry Lamotte Oils/Sana...	Almond oil, refined	130145	1
Lamotte	Almond oil, refined	130145	7
Lamotte/Anzag	Almond oil, refined	130145	1
Lamotte/Sanacorp	Almond oil, refined	130145	1
Henry Lamotte Oils GmbH...	Almond oil, refined	130145	1
AHD	Almond oil, refined	130145	1
Noweda / Henry Lamotte	Almond oil, refined	130145	1
Erika apotheke	Almond oil, refined	130415	1
Eigenherstellung	Almond oil, refined	21112HML01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Almond oil, refined	3236406	1
Caelo	Almond oil, refined	21112PML01	1
Taoasis	Almond oil, refined	3548810	1
Primavera	Almond oil, refined	3808131	2
Caelo	Almond oil, refined	4061202	1
Eigenherstellung	Almond oil, refined	40613HML01	1
Eigenherstellung	Almond oil, refined	40613PML01	1
Eigenherstellung	Almond oil, refined	40613PML02	1
Primavera	Almond oil, refined	4811131	1
Eigenherstellung	Almond oil, refined	50612HML01	1
Caelo	Almond oil, refined	5111208	1
Caelo	Almond oil, refined	514Q1950	1
Eigenherstellung	Almond oil, refined	60913HJK01	1
Caelo	Almond oil, refined	6081202	1
Caelo	Almond oil, refined	7091206	1
Eigenherstellung	Almond oil, refined	80213HJK02	1
Caelo	Almond oil, refined	11227406	1
Caelo	Almond oil, refined	914E1950	1
Lamotte/Gehe	Almond oil, refined	1000003242	1
Caelo	Almond oil, refined	1113E-01950	1
Caelo	Almond oil, refined	11227407	1
Caelo	Almond oil, refined	1131303	1
Caelo	Almond oil, refined	11314206	2
Caelo	Almond oil, refined	11314201	1
Noweda	Almond oil, refined	11314209	1
Caelo	Almond oil, refined	11314208	1
Caelo	Almond oil, refined	11314209	3
Caelo	Almond oil, refined	11314210	1
Taoasis	Almond oil, refined	118979-108698	1
Taoasis	Almond oil, refined	119249-109910	1
Phönix	Almond oil, refined	12041801	2
Caelo	Almond oil, refined	12041801	10
Caelo	Almond oil, refined	12041802	1
Caelo	Almond oil, refined	12041803	3
Noweda	Almond oil, refined	12041803	2
Henry Lamotte	Almond oil, refined	12061302	1
Euro OTC	Almond oil, refined	1207014-01	3
Caelo	Almond oil, refined	12080801	3
Noweda	Almond oil, refined	12080802	1
Caelo	Almond oil, refined	12080802	3
Caelo	Almond oil, refined	12080803	7
Noweda	Almond oil, refined	12080803	2
Sanacorp	Almond oil, refined	12080803	1
Caelo	Almond oil, refined	12080806	10
Noweda	Almond oil, refined	12080806	1
Phönix Gotha	Almond oil, refined	12153501	1
Caelo	Almond oil, refined	12153501	7
Caelo	Almond oil, refined	12153504	1
Caelo	Almond oil, refined	12153505	8
phönix	Almond oil, refined	12153505	1
Caelo	Almond oil, refined	12153506	4
Caelo	Almond oil, refined	12357101	17
Noweda	Almond oil, refined	12357101	2
Anzag 09.02.2013	Almond oil, refined	12357101	1
Caelo	Almond oil, refined	12357102	6
Sanacorp	Almond oil, refined	12357102	3
Caelo	Almond oil, refined	12357103	7
Caelo	Almond oil, refined	12357110	1
Fagron	Almond oil, refined	12A18-N09	1
Caelo	Almond oil, refined	12367101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Almond oil, refined	12J26-N05	1
Fagron	Almond oil, refined	12H23-N01	1
J.H.Müller/direkt/17.05...	Almond oil, refined	130103	1
J.H.Müller/direkt/01.08...	Almond oil, refined	130327	2
Euro OTC	Almond oil, refined	1301001-02	2
Caelo	Almond oil, refined	13091308	1
Caelo	Almond oil, refined	13101203	1
Caelo	Almond oil, refined	1312Q-01950	1
Caesar & Loretz GmbH	Almond oil, refined	13178002	1
Caelo	Almond oil, refined	13178003	18
Noweda	Almond oil, refined	13178003	1
Phönix , 03.09.13	Almond oil, refined	13178003	1
Gehe	Almond oil, refined	13178004	1
Caelo	Almond oil, refined	13178009	2
Caelo	Almond oil, refined	13178011	14
Caelo	Almond oil, refined	132363403	1
Caelo	Almond oil, refined	13236401	16
Gehe	Almond oil, refined	13236403	1
Caelo	Almond oil, refined	13236403	8
Caelo	Almond oil, refined	13236404	4
Caelo	Almond oil, refined	13236406	9
Caesar & Loretz GmbH	Almond oil, refined	13236404	1
Sanacorp	Almond oil, refined	13236406	1
Noweda	Almond oil, refined	13236411	2
Caelo	Almond oil, refined	1326403	1
Krieger	Almond oil, refined	13350006	1
Sanacorp	Almond oil, refined	13350007	1
Caesar & Loretz GmbH	Almond oil, refined	13350007	1
Caelo	Almond oil, refined	13H08-N11	1
Fagron	Almond oil, refined	13C25-N20	1
Noweda	Almond oil, refined	13H08-N11	1
Eigenherstellung	Almond oil, refined	140113HML01	1
Caelo	Almond oil, refined	140113PML01	1
Caelo	Almond oil, refined	1404206	1
Bombastus	Almond oil, refined	14139305	1
Caelo	Almond oil, refined	14139301	1
Primavera	Almond oil, refined	15510141	1
Primavera	Almond oil, refined	1551041	2
Caelo	Almond oil, refined	1601311950	1
Caelo	Almond oil, refined	1602211950	1
Eigenherstellung	Almond oil, refined	160414HML01	1
Henry Lamotte Oils GmbH	Almond oil, refined	16092013B	1
Caelo	Almond oil, refined	17021405	1
Caelo	Almond oil, refined	17121303	1
Caelo	Almond oil, refined	2002311950	1
Caelo	Almond oil, refined	20081301	1
Caelo	Almond oil, refined	2013U01950	1
Caelo	Almond oil, refined	20150612-2	1
Caelo	Almond oil, refined	21011304	1
Caelo	Almond oil, refined	21021303	1
Eigenherstellung	Almond oil, refined	210313HJK01	2
Caelo	Almond oil, refined	22071303	1
Caelo	Almond oil, refined	22081205	1
Caelo	Almond oil, refined	22111301	1
Eigenherstellung	Almond oil, refined	240214HJK01	1
Bombastus	Almond oil, refined	276542	3
Caelo	Almond oil, refined	2405211950	1
Caelo	Almond oil, refined	25041303	1
Caelo	Almond oil, refined	2713U1950	1
Gehe 06.05.2013	Almond oil, refined	276611	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Almond oil, refined	277912	2
Bombastus	Almond oil, refined	279230	1
Lamotte / Gehe	Almond oil, refined	3968/V03/14	1
Bombastus	Almond oil, refined	281765	1
Caelo	Almond oil, refined	31121205	1
H. Lamotte	Almond oil, refined	5120901	1
Henry Lamotte	Almond oil, refined	5177001	3
Caelo	Almond oil, refined	4313E1950	1
Lamotte/Anzag 30.06.2012	Almond oil, refined	5162801	2
Lamotte/Anzag 12.06.2012	Almond oil, refined	5177001	1
Noweda	Almond oil, refined	5177001	4
H. Lamotte	Almond oil, refined	5177001	3
Henry Lamotte/Anzag	Almond oil, refined	5177001	1
henry Lamotte	Almond oil, refined	5177001	1
Lamotte/Krieger	Almond oil, refined	5177001	1
ANZAG	Almond oil, refined	5214206	1
Caelo	Almond oil, refined	5214206	3
Lamotte/Sanacorp	Almond oil, refined	5214206	1
Phönix/Lamotte	Almond oil, refined	5214206	3
Lamotte/gehe	Almond oil, refined	5214206	1
Noweda / Henry Lamotte	Almond oil, refined	5214206	3
Lamotte Phönix	Almond oil, refined	5214206	1
Lamotte	Almond oil, refined	5214206	3
Kehr / Lamotte	Almond oil, refined	5214206	1
H.Lamotte/Noweda	Almond oil, refined	5214206	1
Sanacorp	Almond oil, refined	5214206	2
Henry Lamotte	Almond oil, refined	5214206	1
Henry Lamotte Oils GmbH	Almond oil, refined	5214206	1
Anzag	Almond oil, refined	5214206	1
Henry Lamotte Oils	Almond oil, refined	5214206	1
Henry Lamotte Oils GmbH	Almond oil, refined	5269301	4
Noweda	Almond oil, refined	5269301	4
Klenk	Almond oil, refined	5269301	2
Lamotte/Krieger	Almond oil, refined	5269301	2
Henry Lamotte GmbH	Almond oil, refined	5269301	1
HenryLamotte/Phoenix	Almond oil, refined	5269301	2
H.Lamotte/Noweda	Almond oil, refined	5269301	2
Phönix/Lamotte	Almond oil, refined	5269301	3
Lamotte	Almond oil, refined	5269301	3
Phoenix (Henry Lamotte ...	Almond oil, refined	5269301	1
Phoenix	Almond oil, refined	5269301	1
Lamotte/gehe	Almond oil, refined	5274206	1
Lamotte/Phoenix	Almond oil, refined	5415901	2
Phoenix/Lamotte	Almond oil, refined	5415901	1
Stumpf/Henry Lamotte	Almond oil, refined	5415901	2
Lamotte/ Noweda	Almond oil, refined	5415901	1
Henry Lamotte	Almond oil, refined	5369301	1
LAMOTTE/NOWEDA	Almond oil, refined	5415901	4
Sanacorp	Almond oil, refined	5415901	2
Henry Lamotte	Almond oil, refined	5415901	6
lamotte/anzag	Almond oil, refined	5415901	3
Lamotte/Ebert	Almond oil, refined	5415901	1
Henry Lamotte Oils GmbH	Almond oil, refined	5415901	1
H.Lamotte/Noweda	Almond oil, refined	5415901	1
Henry Lamotte	Almond oil, refined	543802	1
Henry Lamotte/ 11.05.13...	Almond oil, refined	5438301	1
Phönix	Almond oil, refined	543801	1
Henry Lamotte	Almond oil, refined	5438301	4
Großhandel	Almond oil, refined	5438301	5
Phönix/Lamotte	Almond oil, refined	5438301	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda / Henry Lamotte	Almond oil, refined	5438301	3
Henry Lamotte OilsGmbH	Almond oil, refined	5438301	1
Caelo	Almond oil, refined	5438301	4
Lamotte	Almond oil, refined	5438301	6
Lamotte/Sanacorp	Almond oil, refined	5438301	1
Henry Lamotte/Phönix	Almond oil, refined	5438301	1
Phönix	Almond oil, refined	5438301	1
Fischar	Almond oil, refined	5438301	2
Lamotte/Noweda	Almond oil, refined	5438301	1
H.Lamotte/Noweda	Almond oil, refined	5438301	2
Anzag	Almond oil, refined	5438301	2
Phönix, 17.06.2013	Almond oil, refined	5438301	1
Lamotte/ Noweda	Almond oil, refined	5438301	1
lamotte	Almond oil, refined	5438301	1
Henry Lamotte Oils GmbH...	Almond oil, refined	5438301	1
Stumpf/Henry Lamotte	Almond oil, refined	5438302	4
Lamotte/GEHE	Almond oil, refined	5438302	3
Großhandel	Almond oil, refined	5438302	1
Henry Lam./ Phönix	Almond oil, refined	5438302	1
Henry Lamotte	Almond oil, refined	5438302	4
Noweda/Lamotte	Almond oil, refined	5438302	3
Caelo	Almond oil, refined	5438302	3
Sanacorp	Almond oil, refined	5438302	1
AHC03.08.13 EK2,36 100ml	Almond oil, refined	5438302	1
H.Lamotte/Noweda	Almond oil, refined	5438302	2
Lamotte	Almond oil, refined	5438302	7
Henry Lamotte Oils, Gehe	Almond oil, refined	5438302	1
HenryLamotteOilsGmbH/Ph...	Almond oil, refined	5438302	2
Henry Lamotte Oils/Gehe	Almond oil, refined	5438302	2
H.Lamotte/Anzag	Almond oil, refined	5438302	1
Gehe - Lamotte	Almond oil, refined	5438302	1
Henry Lamotte Oils GmbH	Almond oil, refined	5438302	1
noweda/Ichtyol	Almond oil, refined	5438302	2
Lamotte/ Phönix	Almond oil, refined	5438302	1
H.Lamotte / Anzag	Almond oil, refined	5438302	1
Kruger/Lamotte/06.01.14	Almond oil, refined	5438302	1
Lamotte / Alliance Heal...	Almond oil, refined	5438302	1
Henry Lamotte Oils GmBh...	Almond oil, refined	5438302	1
Henry Lamotte Oils / Gehe	Almond oil, refined	548302	1
Taoasis	Almond oil, refined	5452701-111406	1
Henry Lamotte Oils GmbH	Almond oil, refined	5493001	1
Stumpf/Henry Lamotte	Almond oil, refined	5493001	4
Taoasis	Almond oil, refined	5452701-111421	1
Henry Lamotte	Almond oil, refined	5493001	1
Henry Lamotte Oils GmbH...	Almond oil, refined	5493001	1
Lamotte	Almond oil, refined	5493001	2
Henry Lamotte Oils GmbH...	Almond oil, refined	5493001	1
Taoasis	Almond oil, refined	5511701-1113571	1
lamotte/Noweda	Almond oil, refined	5615001	1
Phönix	Almond oil, refined	5669601	1
Lamotte	Almond oil, refined	5703001	1
Spinnrad GmbH	Almond oil, refined	713301223	1
Lamotte/gehe	Almond oil, refined	74760	2
Spinnrad GmbH	Almond oil, refined	913301224	1
xxx	Almond oil, refined		1
Primavera 19.10.15	Evening primrose oil, refined	774J232	1
Lamotte / Gehe	Evening primrose oil, refined	124363	4
Lamotte / Noweda	Evening primrose oil, refined	145581	3
lamotte	Evening primrose oil, refined	145581	3
Lamotte//AHD	Evening primrose oil, refined	145581	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte / Gehe	Evening primrose oil, refined	145802	9
Lamotte / Gehe	Evening primrose oil, refined	150827	4
henry lamotte	Evening primrose oil, refined	150827	3
Phönix	Evening primrose oil, refined	150827	1
Primavera	Evening primrose oil, refined	2611142	1
taoasis/Noweda	Evening primrose oil, refined	37282-115267	1
Caelo	Evening primrose oil, refined	5011508	1
Caelo	Evening primrose oil, refined	7011502	1
Caelo	Evening primrose oil, refined	5061505	1
Caelo	Evening primrose oil, refined	7170717	1
allcura/Noweda	Evening primrose oil, refined		1
Lamotte/AHD	Evening primrose oil, refined	10000001044	1
Lamotte/Noweda	Evening primrose oil, refined	1000003443	1
Taoasis	Evening primrose oil, refined	120161-93968	1
Caelo	Evening primrose oil, refined	10101404	1
Caelo	Evening primrose oil, refined	12031604	1
Euro OTC	Evening primrose oil, refined	1210012-01	13
Caelo	Evening primrose oil, refined	1210012-02	1
Caelo	Evening primrose oil, refined	1281315	1
Euro OTC	Evening primrose oil, refined	1210012-02	12
Euro OTC	Evening primrose oil, refined	1304018-01	10
Euro OTC	Evening primrose oil, refined	1304018-02	12
Großhandel	Evening primrose oil, refined	1304018-02	1
Caelo	Evening primrose oil, refined	13121809	12
Caelo	Evening primrose oil, refined	13121816	15
Caelo	Evening primrose oil, refined	13121812	9
Euro OTC	Evening primrose oil, refined	133339	1
phönix	Evening primrose oil, refined	139798	1
Euro OTC	Evening primrose oil, refined	139798	1
Caelo	Evening primrose oil, refined	14002501	1
Caelo	Evening primrose oil, refined	14002503	7
Euro OTC	Evening primrose oil, refined	1401033-01	7
Phoenix	Evening primrose oil, refined	1401033-01	1
Euro OTC	Evening primrose oil, refined	1401033-02	16
Caelo	Evening primrose oil, refined	14010901	1
Caelo	Evening primrose oil, refined	14010903	2
Caelo	Evening primrose oil, refined	14010904	7
Noweda	Evening primrose oil, refined	14010904	1
Gehe	Evening primrose oil, refined	14010904	1
Phoenix 04.12.14	Evening primrose oil, refined	14010905	1
Caelo	Evening primrose oil, refined	14010907	1
Caelo	Evening primrose oil, refined	14010908	18
Gehe	Evening primrose oil, refined	14010908	1
Caelo	Evening primrose oil, refined	14010910	9
Caelo	Evening primrose oil, refined	14010913	10
Caesar & Loretz GmbH	Evening primrose oil, refined	14010910	1
Caelo	Evening primrose oil, refined	14010912	1
Caelo	Evening primrose oil, refined	14010914	2
Caelo	Evening primrose oil, refined	14010915	12
phönix	Evening primrose oil, refined	14010915	1
Caelo	Evening primrose oil, refined	14010918	5
Caelo	Evening primrose oil, refined	14010980	1
Caelo	Evening primrose oil, refined	1401904	1
Caelo	Evening primrose oil, refined	14010920	1
Caelo	Evening primrose oil, refined	1401908	1
J.H. Müller	Evening primrose oil, refined	140419	1
Caelo	Evening primrose oil, refined	1407002-01	2
Euro OTC	Evening primrose oil, refined	1407002-01	13
Euro OTC	Evening primrose oil, refined	140700201	2
Spirig /Noweda	Evening primrose oil, refined	1407002-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Evening primrose oil, refined	140700202	2
Euro OTC	Evening primrose oil, refined	1407002-02	21
EuroTOC / Noweda	Evening primrose oil, refined	140700202	1
Gehe 5.03.15	Evening primrose oil, refined	1407002-02	1
Phoenix	Evening primrose oil, refined	1407002-02	1
Caelo	Evening primrose oil, refined	1407002-02	1
Noweda	Evening primrose oil, refined	1407002-02	1
Caelo	Evening primrose oil, refined	14256810	1
Caelo	Evening primrose oil, refined	14258601	9
Caelo	Evening primrose oil, refined	14258603	11
Caelo	Evening primrose oil, refined	14258604	4
Phoenix	Evening primrose oil, refined	14258604	1
Caelo	Evening primrose oil, refined	14258606	13
Caelo	Evening primrose oil, refined	14258607	4
Caelo	Evening primrose oil, refined	14258608	13
Fagron	Evening primrose oil, refined	14258608	1
Caelo	Evening primrose oil, refined	14258610	19
Caesar & Loretz GmbH	Evening primrose oil, refined	14258611	1
Caelo	Evening primrose oil, refined	14258611	3
Caelo	Evening primrose oil, refined	14258612	3
Caelo	Evening primrose oil, refined	14258616	15
Caelo	Evening primrose oil, refined	14258617	10
Phoenix	Evening primrose oil, refined	14258617	2
Fagron	Evening primrose oil, refined	14C07-801-291800	1
Caelo	Evening primrose oil, refined	1428616	1
Euro OTC	Evening primrose oil, refined	14C07-B01-291800	1
Fagron	Evening primrose oil, refined	14C07-B01-291800	1
Fagron	Evening primrose oil, refined	14H18-B04-267707	1
Fagron	Evening primrose oil, refined	14H18-B04-297707	3
Ph.Eur.8.2	Evening primrose oil, refined	14H18B04297708	1
Caelo	Evening primrose oil, refined	150081403	1
Euro OTC	Evening primrose oil, refined	150303101	1
Euro OTC	Evening primrose oil, refined	1503031-01	18
Caelo	Evening primrose oil, refined	150303101	1
Caelo	Evening primrose oil, refined	1503031-01	1
	Evening primrose oil, refined	150303101	1
	Evening primrose oil, refined	1503031-01	1
Euro OTC	Evening primrose oil, refined	1506001-01	20
Fagron	Evening primrose oil, refined	1506001-01	1
EuRho	Evening primrose oil, refined	1506001-01	1
Caelo	Evening primrose oil, refined	1506001-01	1
Rondell Apotheke	Evening primrose oil, refined	1506001-01	1
Caelo	Evening primrose oil, refined	150601-01	1
Caelo	Evening primrose oil, refined	15081401	3
Caelo	Evening primrose oil, refined	15081402	3
Caelo	Evening primrose oil, refined	15081403	6
Caelo	Evening primrose oil, refined	15081404	11
Caelo	Evening primrose oil, refined	15081405	18
Caelo	Evening primrose oil, refined	15081407	11
Euro OTC	Evening primrose oil, refined	1509005-01	18
EuRho	Evening primrose oil, refined	1509005-01	2
Caelo	Evening primrose oil, refined	1512000501	1
Euro OTC	Evening primrose oil, refined	1512005-01	21
Euro/Sanacorp	Evening primrose oil, refined	1512005-01	1
Euro OTC	Evening primrose oil, refined	1512005-02	16
Euro/Phoenix	Evening primrose oil, refined	1512005-02	2
Caelo	Evening primrose oil, refined	15121405	1
Caelo	Evening primrose oil, refined	15171101	4
Caelo	Evening primrose oil, refined	15141102	1
Caelo	Evening primrose oil, refined	15171102	19

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Evening primrose oil, refined	15171103	7
Caelo	Evening primrose oil, refined	15171104	15
Caesar & Loretz GmbH	Evening primrose oil, refined	15171104	1
Caelo	Evening primrose oil, refined	15171105	9
Bombastus	Evening primrose oil, refined	15171108	1
Caelo	Evening primrose oil, refined	15171109	9
Caesar & Loretz GmbH	Evening primrose oil, refined	15171109	1
Noweda	Evening primrose oil, refined	15171109	2
Sanacorp	Evening primrose oil, refined	15171109	1
Caelo	Evening primrose oil, refined	15171110	26
Caelo	Evening primrose oil, refined	15171111	12
Caelo	Evening primrose oil, refined	15171113	1
Caelo	Evening primrose oil, refined	15417801	30
Caesar & Loretz GmbH	Evening primrose oil, refined	15417801	1
Caelo	Evening primrose oil, refined	15417802	33
Wolf-Apotheke	Evening primrose oil, refined	15417802	2
Phoenix	Evening primrose oil, refined	15417802	1
Phoenix	Evening primrose oil, refined	15417803	1
Caelo	Evening primrose oil, refined	154178020	1
Caelo	Evening primrose oil, refined	15417803	13
Fagron	Evening primrose oil, refined	15417803	3
Caelo	Evening primrose oil, refined	15417805	12
Caelo	Evening primrose oil, refined	15417806	10
Alte Apotheke	Evening primrose oil, refined	15417808	1
Caelo	Evening primrose oil, refined	15417807	1
Caelo	Evening primrose oil, refined	15417808	6
Caelo	Evening primrose oil, refined	15417809	5
Caelo	Evening primrose oil, refined	154178806	1
Caelo	Evening primrose oil, refined	1541801	1
Fagron	Evening primrose oil, refined	15C02-B99-305094	1
Euro OTC	Evening primrose oil, refined	1603048-01	19
Euroo/Gehe	Evening primrose oil, refined	1603048-01	1
Euro OTC	Evening primrose oil, refined	160304801	2
Caelo	Evening primrose oil, refined	1603048-01	1
Phönix	Evening primrose oil, refined	1603048-01	1
Euro OTC	Evening primrose oil, refined	1605027-01	12
GEHE	Evening primrose oil, refined	1605027-01	1
Caelo	Evening primrose oil, refined	160502701	1
Henry Lamotte	Evening primrose oil, refined	160502701	1
Euro OTC	Evening primrose oil, refined	1608008-01	12
Euro OTC	Evening primrose oil, refined	160800801	1
Caelo	Evening primrose oil, refined	16281301	4
Caelo	Evening primrose oil, refined	16281302	7
Caelo	Evening primrose oil, refined	16281303	2
Caelo	Evening primrose oil, refined	16281305	3
Caelo	Evening primrose oil, refined	16281306	6
Caelo	Evening primrose oil, refined	16281307	7
Caelo	Evening primrose oil, refined	16281308	6
Caelo	Evening primrose oil, refined		3
Caelo	Evening primrose oil, refined	16281309	21
Caelo	Evening primrose oil, refined	16281310	1
Caelo	Evening primrose oil, refined	16281313	4
Caelo	Evening primrose oil, refined	16281318	2
Caelo	Evening primrose oil, refined	162891302	1
Finteler Apotheke	Evening primrose oil, refined	1703020-01	1
Caelo	Evening primrose oil, refined	16281319	1
Euro OTC	Evening primrose oil, refined	1703020-01	13
Caelo	Evening primrose oil, refined	1703020-01	2
Euro OTC	Evening primrose oil, refined	1704018-01	10
Euro OTC	Evening primrose oil, refined	17040118-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Evening primrose oil, refined	1704018-01	1
J.H.Müller	Evening primrose oil, refined	170832	1
Caelo	Evening primrose oil, refined	1710011-01	1
Euro OTC	Evening primrose oil, refined	171001101	1
Caelo	Evening primrose oil, refined	17255801	3
J.H.Müller	Evening primrose oil, refined	171220	1
Caelo	Evening primrose oil, refined	17255802	6
Caelo	Evening primrose oil, refined	17255803	2
Caelo	Evening primrose oil, refined	18009101	5
Caelo	Evening primrose oil, refined	18009102	2
Caelo	Evening primrose oil, refined	18009103	3
Caelo	Evening primrose oil, refined	18009105	2
Caelo	Evening primrose oil, refined	18009107	1
J.H.Müller	Evening primrose oil, refined	180234	1
Caelo	Evening primrose oil, refined	30031501	1
Euro OTC	Evening primrose oil, refined	1813M-03075	1
Taoasis	Evening primrose oil, refined	3788-116147	1
Taoasis	Evening primrose oil, refined	3843-121697	2
Audor Pharma	Evening primrose oil, refined	392-133N	2
Caelo	Evening primrose oil, refined	458528	1
Euro OTC	Evening primrose oil, refined	700775	1
Weltecke GmbH	Evening primrose oil, refined	80011115	1
Audor Pharma	Evening primrose oil, refined	AP140319	2
Audor Pharma	Evening primrose oil, refined	AP140319N	2
Audor Pharma	Evening primrose oil, refined	AP294-118N	10
Audor Pharma	Evening primrose oil, refined	AP392-133N	5
Audor Pharma	Evening primrose oil, refined	AP562-124N	8
Audor Pharma	Evening primrose oil, refined	AP788-128N	9
Audor Pharma	Evening primrose oil, refined	AP141023N	1
Audor Pharma	Evening primrose oil, refined	AP273-116N	5
Audor Pharma	Evening primrose oil, refined	AP396-120N	5
Allianz Healthcare 01.0...	Evening primrose oil, refined	AP562-124N	2
Audor Pharma	Evening primrose oil, refined	AP-788-128N	1
Bergland/Phönix	Evening primrose oil, refined	H302605	1
Bergland/Phönix	Evening primrose oil, refined	K501908	1
Bergland/Anzag	Evening primrose oil, refined	H500690	1
Caelo	Evening primrose oil, refined	122966	3
Allcura	Evening primrose oil, refined	L-387	1
Bergland	Evening primrose oil, refined	K6013776	1
Klenk	Evening primrose oil, refined	122966	4
Lamotte / Gehe	Evening primrose oil, refined	123121	3
Henry Lamotte	Evening primrose oil, refined	123121	1
Lamotte	Evening primrose oil, refined	123121	1
Lamotte / Gehe	Evening primrose oil, refined	123386	7
Lamotte / Gehe	Evening primrose oil, refined	145553	5
Anzag 30.05.2012	Evening primrose oil, refined	144906	1
Lamotte	Evening primrose oil, refined	145553	1
Lamotte / Gehe	Evening primrose oil, refined	145675	4
Anzag/Henry Lamotte	Evening primrose oil, refined	145675	1
Caelo	Evening primrose oil, refined	145691	14
Stumpf/Henry Lamotte	Evening primrose oil, refined	145691	1
Fagron	Evening primrose oil, refined	145691	1
Caelo	Evening primrose oil, refined	2011409	1
Caelo	Evening primrose oil, refined	4061305	1
Euro OTC	Evening primrose oil, refined	80616-11	1
Caelo	Evening primrose oil, refined	11019201	1
Euro OTC	Evening primrose oil, refined	1109004-01	2
Gehe	Evening primrose oil, refined	1109004-02	1
Caesar & Loretz GmbH 40...	Evening primrose oil, refined	11150007	1
Caelo	Evening primrose oil, refined	12032906	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Evening primrose oil, refined	12032908	1
Caelo	Evening primrose oil, refined	12032908	2
Caelo	Evening primrose oil, refined	12032901	1
Noweda	Evening primrose oil, refined	12032910	1
Caelo	Evening primrose oil, refined	12032912	2
Anzag ek:13.04EUR 08.08.12	Evening primrose oil, refined	1205023-01	2
Caelo	Evening primrose oil, refined	12032913	2
Caesar & Loretz GmbH, 4...	Evening primrose oil, refined	12032913	1
Euro OTC	Evening primrose oil, refined	1205023-01	5
Phönix 17.08.2012	Evening primrose oil, refined	1205023-01	2
Noweda	Evening primrose oil, refined	1205023-01	5
Caelo	Evening primrose oil, refined	1205023-02	2
Caelo	Evening primrose oil, refined	120502302	1
Gehe	Evening primrose oil, refined	1205023-02	1
Fagron	Evening primrose oil, refined	1205023-02	1
Euro OTC	Evening primrose oil, refined	1205023-02	7
Anzag 30.01.13 ek:16,40EUR	Evening primrose oil, refined	1205023-02	1
Großhandel	Evening primrose oil, refined	1205023-02	1
Euro OTC	Evening primrose oil, refined	1205023-1	1
Caelo	Evening primrose oil, refined	120543078-13	1
Caelo	Evening primrose oil, refined	1210012-01	2
Jenne	Evening primrose oil, refined	1210012-02	1
EuRho	Evening primrose oil, refined	1210012-02	3
phönix	Evening primrose oil, refined	1210012-02	1
Fiebig 12.07.13 EK15,38...	Evening primrose oil, refined	1210012-02	1
Gehe	Evening primrose oil, refined	121112906	1
Caelo	Evening primrose oil, refined	12112901	2
noveda	Evening primrose oil, refined	12112901	1
ANZAG 27.12.2012	Evening primrose oil, refined	12112901	1
Caelo	Evening primrose oil, refined	12112902	2
Caelo	Evening primrose oil, refined	12112903	1
Noweda	Evening primrose oil, refined	12112905	2
Caelo	Evening primrose oil, refined	12112905	4
Klenk	Evening primrose oil, refined	12112905	1
Kehr	Evening primrose oil, refined	12112906	1
Caelo	Evening primrose oil, refined	12112906	2
Euro OTC	Evening primrose oil, refined	12112906	1
Gehe	Evening primrose oil, refined	12112906	1
Caelo	Evening primrose oil, refined	12112908	5
Caelo	Evening primrose oil, refined	12112910	1
Caelo	Evening primrose oil, refined	12112911	4
Caesar & Loretz GmbH, 4...	Evening primrose oil, refined	12112911	1
Caelo	Evening primrose oil, refined	12112912	1
Ichthyl Phönix	Evening primrose oil, refined	12112912	1
Caelo	Evening primrose oil, refined	12112914	1
Caelo	Evening primrose oil, refined	12112915	3
Hustadt Apotheke	Evening primrose oil, refined	1211905	1
Caelo	Evening primrose oil, refined	12325101	8
Fagron	Evening primrose oil, refined	12325101	2
Noweda	Evening primrose oil, refined	12325101	1
Caesar & Loretz GmbH, 4...	Evening primrose oil, refined	12325101	1
Caelo	Evening primrose oil, refined	12325103	5
Caelo	Evening primrose oil, refined	12325105	7
Caelo	Evening primrose oil, refined	12325106	16
Kehr	Evening primrose oil, refined	12325106	1
Caesar & Loretz GmbH, 4...	Evening primrose oil, refined	12325106	1
Gehe	Evening primrose oil, refined	12325106	1
Caelo	Evening primrose oil, refined	12325107	8
Caesar & Loretz GmbH	Evening primrose oil, refined	12325107	1
Caelo	Evening primrose oil, refined	12325109	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Evening primrose oil, refined	12325110	2
Gehe	Evening primrose oil, refined	12325110	1
Caelo	Evening primrose oil, refined	12325111	11
Sanacorp	Evening primrose oil, refined	12325111	1
Caelo	Evening primrose oil, refined	12325113	7
Caesar & Loretz GmbH, 4. . .	Evening primrose oil, refined	12325113	1
Gehe	Evening primrose oil, refined	12325113	1
Caelo	Evening primrose oil, refined	12325114	4
Sanacorp	Evening primrose oil, refined	12325114	1
Bombastus	Evening primrose oil, refined	12325114	1
Caelo	Evening primrose oil, refined	12325115	3
Caelo	Evening primrose oil, refined	12325116	4
Caelo	Evening primrose oil, refined	1235101	1
Caelo	Evening primrose oil, refined	1235106	1
Caelo	Evening primrose oil, refined	1235113	1
Caelo	Evening primrose oil, refined	1232516	1
Caelo	Evening primrose oil, refined	12355115	1
Fagron	Evening primrose oil, refined	12e23-N06	1
Fagron	Evening primrose oil, refined	12J08-N04	1
J.H.Müller/direkt/17.05. . .	Evening primrose oil, refined	130127	1
Euro/Phönix	Evening primrose oil, refined	1304018-01	2
Euro OTC	Evening primrose oil, refined	130401801	1
Großhandel	Evening primrose oil, refined	1304018-01	1
Phönix 27.01.2014	Evening primrose oil, refined	1304018-01	1
Euro/ Phönix	Evening primrose oil, refined	1304018-02	1
Caelo	Evening primrose oil, refined	131211807	1
Anzag	Evening primrose oil, refined	130401812	1
Caelo	Evening primrose oil, refined	13121804	1
Caelo	Evening primrose oil, refined	13121806	2
Bombastus	Evening primrose oil, refined	13121806	1
Caelo	Evening primrose oil, refined	13121807	9
Noweda	Evening primrose oil, refined	13121807	1
Caelo	Evening primrose oil, refined	13121808	17
Phönix 18.11.2013	Evening primrose oil, refined	13121808	1
Sanacorp	Evening primrose oil, refined	13121809	2
Caesar & Loretz	Evening primrose oil, refined	13121809	1
Caesar & Loretz GmbH	Evening primrose oil, refined	13121812	1
Caelo	Evening primrose oil, refined	13121813	4
Caelo	Evening primrose oil, refined	13121814	7
celo/phönix	Evening primrose oil, refined	13121814	1
Sanacorp	Evening primrose oil, refined	13121814	1
Caeol/Phönix	Evening primrose oil, refined	13121814	1
Cealo / Noweda	Evening primrose oil, refined	13121814	1
Caelo	Evening primrose oil, refined	13121818	4
Anzag	Evening primrose oil, refined	13121818	1
Herbasinica	Evening primrose oil, refined	1312809	2
Fargron	Evening primrose oil, refined	13E06-N09	1
Fagron	Evening primrose oil, refined	13E06-N09	3
Caelo	Evening primrose oil, refined	1312812	1
Cae	Evening primrose oil, refined	14002501	1
Caelo	Evening primrose oil, refined	14002502	2
müller	Evening primrose oil, refined	1401332X	1
Lamotte	Evening primrose oil, refined	145675	1
Euro OTC	Evening primrose oil, refined	1506001	1
dr.r.hahn	Evening primrose oil, refined	1509005-01	1
Caelo	Evening primrose oil, refined	160213003	1
Caelo	Evening primrose oil, refined	20140612-18	1
Caelo	Evening primrose oil, refined	3112M-03075	1
PG-Naturpharma/Sanacorp	Evening primrose oil, refined	711-098	1
weltecke/Noweda	Evening primrose oil, refined	80010212	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Audor Pharma	Evening primrose oil, refined	869-100N	1
Audor Pharma	Evening primrose oil, refined	AP191-088N	2
Audor Pharma	Evening primrose oil, refined	AP920-114N	9
Audor Pharma	Evening primrose oil, refined	AP869-100/1N	1
Audor Pharma	Evening primrose oil, refined	AP094-115N	1
Audor Pharma	Evening primrose oil, refined	AP191-088/1N	1
Audor Pharma	Evening primrose oil, refined	AP366-108N	3
Audor Pharma	Evening primrose oil, refined	AP676-078N	1
Audor Pharma	Evening primrose oil, refined	AP869-100N	1
Andor	Evening primrose oil, refined	AP869-100N	1
Audor Pharma	Evening primrose oil, refined	AP920-114/1N	1
Audor Pharma	Evening primrose oil, refined	Ap869-100N	1
Phoenix	Evening primrose oil, refined	AP869100N	1
Andor	Evening primrose oil, refined	AP920-114N	1
Caelo	Jojoba oil, native	16070609	1
allcura	Evening primrose oil, refined	L2641	1
Caelo	Jojoba oil, native	17073310	2
Lamotte/Fiebig	Jojoba oil, native	8002391001	1
Caelo	Flaxseed oil, native	14306806	1
Caelo	Flaxseed oil, native	15175604	1
Caelo	Flaxseed oil, native	17236401	1
Caelo	Flaxseed oil, native	17236408	1
Caelo	Flaxseed oil, native	263894	1
Euro OTC	Flaxseed oil, native	13058506	1
Caelo	Flaxseed oil, native	13178210	2
statfold	Flaxseed oil, native	8001-26-1	1
Caelo	Olive oil, native	1081402	1
Henry Lamotte Oils GmbH	Olive oil, native	124392	1
Bombastus	Olive oil, native	6861571	2
Caelo	Olive oil, native	5011503	1
Caelo	Olive oil, native	2061503	1
Henry Lamotte/ Noweda	Olive oil, native	124392	1
lamotte/noweda	Olive oil, native	1000004908	1
Fagron	Olive oil, native	129189	1
Fagron	Olive oil, native	130190	2
Caelo	Olive oil, native	13051402	1
PHOENIX 03.03.2014	Olive oil, native	1305014-01	1
Fagron	Olive oil, native	133026	2
Caelo	Olive oil, native	13147004	19
Caelo	Olive oil, native	13350102	18
Caesar & Loretz GmbH Hi...	Olive oil, native	13350102	1
Caelo	Olive oil, native	13350104	7
Fagron	Olive oil, native	13J25-B04	1
Fagron	Olive oil, native	13J25-B04-289441	2
Caelo	Olive oil, native	14024406	5
Caelo	Olive oil, native	14024407	11
Caelo	Olive oil, native	14024408	14
Caelo	Olive oil, native	14024409	8
Euro OTC	Olive oil, native	1406009-01	1
Caelo	Olive oil, native	1406914	1
Caelo	Olive oil, native	14070802	1
Caelo	Olive oil, native	14071504	1
Caelo	Olive oil, native	14118603	20
Phoenix	Olive oil, native	14118603	1
Caelo	Olive oil, native	14118604	12
Bombastus	Olive oil, native	14118604	1
Caelo	Olive oil, native	14118606	7
Phönix	Olive oil, native	14118606	1
Caelo	Olive oil, native	14182301	4
Caelo	Olive oil, native	14182302	16

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda 21.02.2015	Olive oil, native	14182302	1
Caelo	Olive oil, native	14182303	21
Caelo	Olive oil, native	14182304	7
Euro OTC	Olive oil, native	14182303	1
Caelo	Olive oil, native	14306902	1
Caelo	Olive oil, native	14306903	4
Caelo	Olive oil, native	14306904	7
Bombastus	Olive oil, native	14306904	1
Caesar & Loretz GmbH	Olive oil, native	14306904	1
sanacorp	Olive oil, native	14306904	1
Caelo	Olive oil, native	14306905	13
Caelo	Olive oil, native	14306906	3
Caelo	Olive oil, native	14306915	7
Caelo	Olive oil, native	14306914	1
Caelo	Olive oil, native	14306916	11
Caelo	Olive oil, native		3
Fagron	Olive oil, native	14B12-B06	1
Sanacorp 21.01.15 EK 10...	Olive oil, native	14G08-B02-300073	1
Caelo	Olive oil, native	15112703	3
Fagron	Olive oil, native	14G08-B02-300074	3
Caelo	Olive oil, native	15112705	7
Caelo	Olive oil, native	15112704	2
sanacorp	Olive oil, native	15112706	1
Caelo	Olive oil, native	15112706	13
Caelo	Olive oil, native	15112707	4
Caelo	Olive oil, native	15216307	10
Caelo	Olive oil, native	15216308	11
Phönix Gotha	Olive oil, native	15216308	1
Caelo	Olive oil, native	15216309	7
Caesar & Loretz GmbH Hi...	Olive oil, native	15216309	1
Caelo	Olive oil, native	15216310	4
Caelo	Olive oil, native	15311401	6
Caelo	Olive oil, native	15311402	5
Caelo	Olive oil, native	15311403	12
Caelo	Olive oil, native	15311404	10
Caelo	Olive oil, native	15311408	6
Caelo	Olive oil, native	15399702	12
Bombastus	Olive oil, native	15399702	1
Noweda	Olive oil, native	15399702	1
Caelo	Olive oil, native	15399703	13
Caelo	Olive oil, native	15399704	11
Caesar&lorenz,Phoenix,1...	Olive oil, native	15399703	1
Fagron	Olive oil, native	15A27-B03-304709	5
Fagron	Olive oil, native	15A27-B03	1
Caelo	Olive oil, native	16035308	6
Caelo	Olive oil, native	15399706	1
Caelo	Olive oil, native	16035309	4
Caelo	Olive oil, native	16035310	13
Phoenix	Olive oil, native	16035310	1
Caelo	Olive oil, native	16035311	6
Phoenix	Olive oil, native	16035311	1
Caelo	Olive oil, native	16144605	3
Caelo	Olive oil, native	16212801	7
Gehe	Olive oil, native	16212801	1
Caelo	Olive oil, native	16212802	2
Caelo	Olive oil, native	16212803	4
Caelo	Olive oil, native	16212805	1
Caelo	Olive oil, native	16212806	17
Caelo	Olive oil, native	16278806	1
Henry Lamotte	Olive oil, native	16327402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PHÖNIX	Olive oil, native	16327402	1
Caesar&LoretzGmbH,07,03...	Olive oil, native	16327402	1
Caelo	Olive oil, native	16327402	12
Phönix	Olive oil, native	16327402	1
Allianz Healthcare	Olive oil, native	16327402	1
Caesar 6 Loretz GmbH	Olive oil, native	16327402	1
Caelo	Olive oil, native	16327403	6
Caelo	Olive oil, native	16327406	7
Phönix	Olive oil, native	16327406	1
Fagron	Olive oil, native	16K24-01-329002	1
Caelo	Olive oil, native	16K24-B01329002	1
Phönix	Olive oil, native	17038702	1
Caelo	Olive oil, native	17020700/16212806	1
Caelo	Olive oil, native	17038702	13
Caesar&Loretz	Olive oil, native	17038702	2
Caelo	Olive oil, native	17038703	7
Caesar & Loretz GmbH,20...	Olive oil, native	17121602	1
Caelo	Olive oil, native	17121602	5
Caelo	Olive oil, native	17121603	1
Caelo	Olive oil, native	17121604	8
Caelo	Olive oil, native	17121606	5
Caelo	Olive oil, native	17121900/17121610	1
Caelo	Olive oil, native	17236501	3
Phönix	Olive oil, native	17121610	1
Sanacorp	Olive oil, native	17236501	1
Caelo	Olive oil, native	17236505	5
Caelo	Olive oil, native	17236502	1
Caelo	Olive oil, native	17361101	5
Caelo	Olive oil, native	17236506	1
GEHE	Olive oil, native	18002118	1
Bombastus	Olive oil, native	18002119	3
Caelo	Olive oil, native	20061502	1
Caelo	Olive oil, native	24031502	1
Caelo	Olive oil, native	24081507	1
Bombastus	Olive oil, native	279608	5
Bombastus	Olive oil, native	282963	7
gehe	Olive oil, native	282963	2
gehe	Olive oil, native	284144	2
Bombastus	Olive oil, native	284144	8
gehe	Olive oil, native	285942	3
Bombastus	Olive oil, native	285942	10
Bombastus	Olive oil, native	286533	8
Bombastus	Olive oil, native	287259	7
Gehe	Olive oil, native	287259	1
Brücken-Apotheke	Olive oil, native	288752	2
Bombastus	Olive oil, native	288752	5
Caelo	Olive oil, native	288752	2
Bombastus	Olive oil, native	290403	10
Caelo	Olive oil, native	290924	1
Bombastus	Olive oil, native	289095	5
Bombastus	Olive oil, native	290924	17
Bombastus	Olive oil, native	291374	7
Beiersdorf AG/Geilenkir...	Olive oil, native	291374	1
Caelo	Olive oil, native	2915A-02000	1
Bombastus	Olive oil, native	292369	16
Caelo	Olive oil, native	292369	2
Bombastus	Olive oil, native	294034	13
Caelo	Olive oil, native	294034	1
Bombastus	Olive oil, native	295282	4
Bombastus	Olive oil, native	295283	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Olive oil, native	296259	2
Bombastus	Olive oil, native	296260	12
Bombastus	Olive oil, native	297564	2
Bombastus	Olive oil, native	298008	7
Bombastus	Olive oil, native	297227	5
Bombastus	Olive oil, native	297565	3
Bombastus	Olive oil, native	298007	4
Bombastus	Olive oil, native	300219	5
Bombastus	Olive oil, native	300220	10
Bombastus	Olive oil, native	301100	16
Bombastus	Olive oil, native	301099	6
Bombastus	Olive oil, native	302522	3
Bombastus	Olive oil, native	3011000	2
Bombastus	Olive oil, native	302523	9
Bombastus	Olive oil, native	304112	1
Bombastus	Olive oil, native	304113	5
Bombastus	Olive oil, native	305002	7
GEHE	Olive oil, native	305002	1
Bombastus	Olive oil, native	305003	9
Bombastus	Olive oil, native	306949	5
Bombastus	Olive oil, native	306950	7
Caelo	Olive oil, native	306950	1
Bombastus	Olive oil, native	307103	3
Caelo	Olive oil, native	307103	3
Bombastus	Olive oil, native	308230	13
Bombastus	Olive oil, native	308229	3
GEHE	Olive oil, native	308229	1
Bombastus	Olive oil, native	310868	6
Caelo	Olive oil, native	3916I-02000	1
Raben-Apotheke Mildenau	Olive oil, native	5623601	1
Gehe - Lamotte	Olive oil, native	5470202	1
Hedinger	Olive oil, native	5623601	1
Lamotte / Alliance Heal...	Olive oil, native	5671601	1
Henry Lamotte	Olive oil, native	5671601	1
Henry Lamotte	Olive oil, native	5671901	3
H. Lamotte/Noweda	Olive oil, native	5671901	1
Lamotte Gmbh	Olive oil, native	5671901	2
Lamotte/Sanacorp	Olive oil, native	5671901	1
Caelo	Olive oil, native	5734701	1
Henry Lamotte	Olive oil, native	5734701	1
Noweda/Lamotte	Olive oil, native	5813601	1
Henry Lamotte	Olive oil, native	5813601	1
Lamotte/Anzag	Olive oil, native	5813601	1
Caelo	Olive oil, native	5891202	1
Noweda-Henry	Olive oil, native	5891202	1
Raben-Apotheke Mildenau	Olive oil, native	5891202	3
Henry Lamotte/Gehe	Olive oil, native	5903400003	1
Caelo	Olive oil, native	60081503	2
Henry Lamotte	Olive oil, native	8000652001	1
Lamotte/Sanacorp	Olive oil, native	8000652001	1
Henry Lamotte/Kehr	Olive oil, native	8000652001	1
Henry Lamotte/Noweda	Olive oil, native	8000821001	1
Henry Lamotte Oils GmbH	Olive oil, native	8001560001	1
Henry Lamotte	Olive oil, native	8000821001	1
Henry Lamotte/Phoenix	Olive oil, native	8001431001	1
Henry Lamotte	Olive oil, native	8001560001	1
Henry Lamotte/Phönix	Olive oil, native	8001648001	1
Caelo	Olive oil, native	8001854001	2
lamotte/noweda	Olive oil, native	8001854001	1
Henry Lamotte/Phoenix	Olive oil, native	8001948002	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte Gmbh	Olive oil, native	8001948002	1
Henry Lamotte/Phönix	Olive oil, native	8001948002	1
Gehe-Henry	Olive oil, native	8002026001	1
Gehe	Olive oil, native	8002366001	1
lamotte/noweda	Olive oil, native	8002827001	1
Henry Lamotte/Noweda	Olive oil, native	8002827001	1
Henry Lamotte Oils GmbH	Olive oil, native	8003420001	1
Lamotte/Kehr	Olive oil, native	8003580002	1
Henry Lamotte/Sanacorp	Olive oil, native	8003801001	1
Henry Lamotte Oils GmbH...	Olive oil, native	8003801001	1
Bombastus	Olive oil, native		1
Wolfvital	Olive oil, native	W48B041	1
phönix	Olive oil, native		1
Noweda	Olive oil, native	272270	2
Bombastus	Olive oil, native	273871	5
Bombastus	Olive oil, native	274236	10
Caelo	Olive oil, native	114Q-02000	1
Lamotte	Olive oil, native	123773	2
Henry Lamotte/Anzag	Olive oil, native	123773	1
Henry Lamotte / Noweda	Olive oil, native	123773	1
Caelo	Olive oil, native	3071302	1
Caelo	Olive oil, native	11376905	3
Caelo	Olive oil, native	11376906	4
Caelo	Olive oil, native	12077701	17
Caelo	Olive oil, native	12077702	4
fargon	Olive oil, native	1212022-01	2
Caelo	Olive oil, native	12141305	3
Caelo	Olive oil, native	12141306	11
Caelo	Olive oil, native	12141309	6
Caelo	Olive oil, native	12274104	7
Caelo	Olive oil, native	12274105	1
Caelo	Olive oil, native	12274106	12
Alliance Healthcare	Olive oil, native	12274108	1
Bombastus	Olive oil, native	12274106	1
Euro OTC	Olive oil, native	12274109	3
Fagron	Olive oil, native		1
Caelo	Olive oil, native	13010904	24
celo/Jenne	Olive oil, native	13010904	1
phönix	Olive oil, native	13010904	1
Caelo	Olive oil, native	13010906	11
Gehe 13.06.2013	Olive oil, native	13010906	2
Caelo	Olive oil, native	13010907	1
Caelo	Olive oil, native	130109062x	1
Caelo	Olive oil, native	13010911	19
EuRho	Olive oil, native	13010911	1
Caelo	Olive oil, native	13147001	10
Gehe	Olive oil, native	13010911	1
Caelo	Olive oil, native	13147002	9
sanacorp	Olive oil, native	13147004	1
Caelo	Olive oil, native	13224701	6
Gehe 04.11.2013	Olive oil, native	13224701	1
Caelo	Olive oil, native	13224702	4
Caelo	Olive oil, native	13224703	20
Caelo	Olive oil, native	13224704	16
Caelo	Olive oil, native	13350101	4
Noweda	Olive oil, native	13350102	1
Caelo	Olive oil, native	13350105	3
Fagron	Olive oil, native	13C08-N06	8
Sanacorp	Olive oil, native	13C08-N06	1
Fagron	Olive oil, native	13G18-B08	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Olive oil, native	13G18-B08-286424	3
Caelo	Olive oil, native	16011404	1
Caelo	Olive oil, native	25021404	1
Phönix	Olive oil, native	20032	1
Caelo	Olive oil, native	16121306	1
Bombastus	Olive oil, native	275757	3
Phönix Gotha	Olive oil, native	277049	1
Noweda	Olive oil, native	277050	2
Gehe	Olive oil, native	277050	2
Bombastus	Olive oil, native	278583	9
Bombastus	Olive oil, native	277050	9
Bombastus	Olive oil, native	279609	7
Bombastus	Olive oil, native	281209	6
Bombastus	Olive oil, native	281898	6
Caelo	Olive oil, native	28121202	1
Gehe	Olive oil, native	282963	2
Lamotte/Sanacorp	Olive oil, native	5271504	2
Caelo	Olive oil, native	31081206	1
Caelo	Olive oil, native	30071305	1
Lamotte/ Anzag	Olive oil, native	5271504	1
H. Lamotte	Olive oil, native	5271504	1
Lamotte Phönix	Olive oil, native	5271504	1
Caelo	Olive oil, native	5271504	2
lamotte/Gehe	Olive oil, native	5271504	1
Gehe - Lamotte	Olive oil, native	5271504	1
Noweda/Anzag	Olive oil, native	5272401	1
Henry Lamotte	Olive oil, native	5272401	3
Lamotte	Olive oil, native	5272401	3
Henry Lamotte, Noweda	Olive oil, native	5372701	2
Lamotte	Olive oil, native	5464901	2
Henry Lamotte	Olive oil, native	5464901	1
Gehe	Olive oil, native	5470202	1
Lamotte/ Anzag	Olive oil, native	5470202	2
Lamotte	Olive oil, native	5470202	3
Phoenix	Olive oil, native	5470202	1
Lamotte / Gehe	Olive oil, native	5470202	1
H.Lamotte/Noweda	Olive oil, native	5470202	1
Caelo	Olive oil, native	90417099	1
Euro OTC	Evening primrose oil, refined	140103302	1
Henry Lamotte/Anzag	Evening primrose oil, refined	124363	2
Hedinger	Evening primrose oil, refined	14010915	5
Caelo	Evening primrose oil, refined	14010917	1
Caelo	Evening primrose oil, refined	14258615	2
Euro OTC	Evening primrose oil, refined	13040-18-01	1
Caelo	Evening primrose oil, refined	14010905	1
VDL;22.05.15;11.47EUR	Almond oil, refined	14215807	1
Caelo	Almond oil, refined	15127303	2
Sanacorp, 4,88EUR, 18.1.17	Almond oil, refined	16327303	1
Sanacorp, 11.03.17, 5,46EUR	Almond oil, refined	16327309	2
Henry Lamotte Oils GmbH...	Almond oil, refined	1714U-01950	1
Caelo	Almond oil, refined	3714U-01950	1
Hanry Lamotte Oils / No...	Almond oil, refined	8000915001	1
Henry LAmotte Oils GmbH	Almond oil, refined	8001497001	6
Henry Lamotte / Noweda	Almond oil, refined	8000915002	2
HENRY LAMOTTE OILS GmbH	Almond oil, refined	8001497001	1
Henry Lamotte / Noweda	Almond oil, refined	8002653001	1
Noweda	Almond oil, refined	11078806	1
Henry Lamotte / Sanacorp	Almond oil, refined	8001497001	2
Lamotte/Krieger	Almond oil, refined	5081301	1
Lamotte/N/31.05.13/PZN:...	Almond oil, refined	5214206	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte/N/18.06.1...	Almond oil, refined	5438301	1
Sanacorp;13.05.15; 2,94EUR	Peanut oil, refined	14355905	1
Noweda	Peanut oil, refined	11B16-N02	2
Caelo	Peanut oil, refined	12203806	3
Henry Lamotte GmbH/Hold...	Peanut oil, refined	5231601	1
Caelo	Avocado oil, refined	12064812	2
Bergland	Borage oil, refined, stabili...	H300590	2
VDL	Codliver oil (type A)	5878701	1
Lamotte/Noweda	Codliver oil (type A)	8000601001	3
Lamotte/ Noweda	Codliver oil (type A)	8000918001	2
Henry Lamotte;VDL; 22,5...	Codliver oil (type A)	8000601001	1
Henry Lamotte / Anzag	Codliver oil (type A)	8000918001	2
Henry Lamotte GmbH/Hold...	Codliver oil (type A)	5228001	2
Anzag	Codliver oil (type A)	5228001	2
Caelo	Jojoba oil, native	14052206	1
Caelo	Jojoba oil, native	2614A-01986	1
Bombastus	Jojoba oil, native	267834	1
Caelo	Jojoba oil, native	11078606	1
Caelo	Walnut oil, refined	10328802	1
henry lamotte oils	Jojoba oil, native	5055401	1
Caelo	Walnut oil, refined	11167103	1
Henry Lamotte Oils / No...	Flaxseed oil, native	1000000878	1
VDL;2.51EUR;22.06.15	Flaxseed oil, native	15116702	1
Lamotte/AHD	Evening primrose oil, refined	145581	2
Henry Lamotte /AHD	Evening primrose oil, refined	1000001044	1
Bombastus	Evening primrose oil, refined	1000003443	3
Lamotte/AHD	Evening primrose oil, refined	1000004385	1
Lamotte/AHD	Evening primrose oil, refined	1000007046	1
Henry Lamotte /AHD	Evening primrose oil, refined	1000005123	1
Taoasis	Evening primrose oil, refined	3723-107968	1
Taoasis	Evening primrose oil, refined	12096-99883	2
Noweda	Evening primrose oil, refined	E1106072	1
Caelo	Olive oil, native	14306913	1
VDL Nachlieferung; 10.7.15...	Olive oil, native	291374	2
Caelo	Olive oil, native	295282	3
Caelo	Olive oil, native	300117/A	1
Lamotte/Spangro	Olive oil, native	5903400003	1
Lamotte/Spangro	Olive oil, native	5891202	1
Lamotte/Spangro	Olive oil, native	8000652001	1
Lamotte / Phönix	Olive oil, native	8001431001	1
Lamotte/Spangro	Olive oil, native	8001648001	1
Bombastus	Olive oil, native	270970	4
Anzag	Olive oil, native	12077701	1
Fiebig	Olive oil, native	12077702	1
Anzag	Olive oil, native	12274104	1
Gehe	Olive oil, native	279609	1
Gehe	Olive oil, native	288752	2
Bombastus	Olive oil, native	oo270970	1
Caelo	Olive oil, native	5272401	4
HL/N/25.02.13 PZN:3649899	Olive oil, native	5283901	1
HENRY LAMOTTE OILS GmbH	Avocado oil, refined	1000001451	1
VDL;11.3.16;2,77EUR	Avocado oil, refined	15268203	1
Sanacorp, 16.02.17, 3,04EUR	Avocado oil, refined	16104702	1
Sanacorp;31.03.17, 5,20EUR	Avocado oil, refined	16185908	1
lamotte	Sesame oil, refined	123488	1
Caelo	Sesame oil, refined	13172906	6
Caelo	Sesame oil, refined	14170004	6
Caelo	Sesame oil, refined	14170005	3
Caelo	Sesame oil, refined	16312505	1
Henry Lamotte	Sesame oil, refined	5744101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sesame oil, refined	15374001	4
Lamotte/Phönix	Sesame oil, refined	58772901	1
henry lamotte	Sesame oil, refined	122972	1
Henry Lamotte Oils	Sesame oil, refined	122972	1
Caelo	Sesame oil, refined	13011105	10
Caelo	Sesame oil, refined	5319201	2
VDL;11.3.16;3,45EUR	Wheatgerm oil, refined	15374101	2
Caelo	Wheatgerm oil, refined	16137403	3
Caelo	Grape seed oil, refined	13174601	2
Caelo	Grape seed oil, refined	13065302	7
Caelo	Grape seed oil, refined	14127906	1
Caelo	Grape seed oil, refined	15144502	1
Caesar & Lorenz	Grape seed oil, refined	13065310	1
Caelo	Oleyl oleate (stabilised)	14129007	9
Caelo	Oleyl oleate (stabilised)	15410207	3
Sanacorp	Oleyl oleate (stabilised)	13378403	1
Sanacorp	Oleyl oleate (stabilised)	12220208	1
Sanacorp	Oleyl oleate (stabilised)	13096402	2
Caelo	Oleyl oleate (stabilised)	13291602	1
Caelo	Oleyl oleate (stabilised)	13096801	1
Caelo	Oleyl oleate (stabilised)	13291605	1
Caelo	Oleyl oleate (stabilised)	4170717	1
Fagron	Oleyl oleate (stabilised)	11I08-N03	1
Noweda	Oleyl oleate (stabilised)	11/08-N03	1
Fagron	Oleyl oleate (stabilised)	11K25-N09	1
Fagron	Oleyl oleate (stabilised)	11L08-N03	1
Noweda	Oleyl oleate (stabilised)	12067704	1
AHD	Oleyl oleate (stabilised)	14129003	1
Caelo	Oleyl oleate (stabilised)	14129003	1
Phoenix	Oleyl oleate (stabilised)	14129007	1
Caelo	Oleyl oleate (stabilised)	14129007(2mal)	1
Caelo	Oleyl oleate (stabilised)	14254603	7
Fagron	Oleyl oleate (stabilised)	14B14-B03-29219	1
Fagron	Oleyl oleate (stabilised)	14B14-B03-29295	1
Cae	Oleyl oleate (stabilised)	15321701	1
Caelo	Oleyl oleate (stabilised)	15321706	4
Caesar&Loretz/Noweda	Oleyl oleate (stabilised)	15321702	1
Caelo	Oleyl oleate (stabilised)	15321702	1
Caelo	Oleyl oleate (stabilised)	15322703	2
Caelo	Oleyl oleate (stabilised)	15410204	4
Caesar&Loretz GmbH	Oleyl oleate (stabilised)	15410204	1
Caelo	Oleyl oleate (stabilised)	16104907	3
Hubertus Apotheke	Oleyl oleate (stabilised)	15410207	1
Noweda	Oleyl oleate (stabilised)	17013008	1
Caelo	Oleyl oleate (stabilised)	17051700/15410207	1
Caelo	Oleyl oleate (stabilised)	17121500/17013008	1
Caelo	Oleyl oleate (stabilised)	17358101	1
Caelo	Oleyl oleate (stabilised)	17358103	1
Caelo	Olive oil, native	14182305	2
Bombastus	Olive oil, native	287259Noweda	1
Bombastus	Olive oil, native	289094	5
Bombastus	Olive oil, native	294035	2
Bombastus	Olive oil, native	4328413	1
Lamotte/Spangro	Olive oil, native	5671901	1
Pharmachem/AHD	Olive oil, native	13324703	1
Bombastus	Olive oil, native	276542	1
Bombastus	Olive oil, native	277049	1
Sanacorp	Olive oil, native	277050	1
Sanacorp	Olive oil, native	279609	1
Caelo	Olive oil, native	279609	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Olive oil, native	283762	3
Caelo	Olive oil, native	12388	1
Henry lamotte/N/PZN:004...	Olive oil, native	5372701	1
Fagron	Olive oil, native	139812	1
Ebert & Jacobi 02.08...	Olive oil, native	15399704	1
Caesar & Loretz GmbH Hi...	Olive oil, native	16212801	1
Caesar & Loretz GmbH Hi...	Olive oil, native	16212806	1
Frago/Noweda	Olive oil, native	16K24-B01	1
Frago/Gehe	Olive oil, native	16K24-B01-329002	2
Fagron	Olive oil, native	16K24-B01-329002	2
Caelo	Olive oil, native	17121607	1
Caelo	Olive oil, native	17121610	1
Caelo	Olive oil, native	172236501	1
Caelo	Olive oil, native	17236508	1
Caelo	Olive oil, native	17361104	1
Nachprüfung	Olive oil, native	292369	2
Caelo	Olive oil, native	297564	1
Allianc. Health. /bo	Olive oil, native	302523	1
Bomastus/Phönix	Olive oil, native	302523	1
Caelo	Olive oil, native	305003	1
Bombastus	Olive oil, native	310866	1
Fagron	Olive oil, native	329002	1
Henry Lamotte	Olive oil, native	5891202	1
Apotheke am Landrain	Olive oil, native		1
Caelo	Olive oil, native	6/AKD15	1
Gehe/henry	Olive oil, native	8001431001	1
Henry Lamotte/Sanacorp	Olive oil, native	8000652001	1
Lamotte	Olive oil, native	8000821001	1
Lamotte	Olive oil, native	5903400003	1
hen	Olive oil, native	8001648001	1
Henry Lamotte/Anzag	Olive oil, native	8001648001	1
Henry Lamotte	Olive oil, native	8001648001	1
Caelo	Olive oil, native	800236001	1
Lamotte/Noweda	Olive oil, native	8002656001	1
Henry Lamotte/Phönix	Olive oil, native	8003420001	2
Gehe-Henry	Olive oil, native	8002366001	1
Lamotte/Noweda	Olive oil, native	8003580002	1
Henry Lamotte Oils GmbH...	Olive oil, native	124392	1
Lamotte / Noweda	Olive oil, native	124392	1
Caelo	Olive oil, native	2021501	1
Caelo	Olive oil, native	2021502	1
Caelo	Olive oil, native	7051507	1
Caelo	Olive oil, native	7031501	1
Caelo	Olive oil, native	8121404	1
Sanacorp	Olive oil, native	1/AKD4	2
Sanacorp 18.03.14 EK 10...	Olive oil, native	13G18-B08-286423	1
Caelo	Olive oil, native	14024410	4
Caesar & Loretz GmbH Hi...	Olive oil, native	14182302	1
Sanacorp	Olive oil, native	14306903	1
Caelo	Olive oil, native	15216305	1
anzag	Olive oil, native	15216308	1
Caol / Gehe	Olive oil, native	15311401	1
Caolo/Phönix	Olive oil, native	15311402	1
Caesar & Loretz GmbH Hi...	Olive oil, native	15311404	1
Caelo	Olive oil, native	15363302	1
Caesar&Loretz GmbH, AHC...	Olive oil, native	17038702	1
Caelo	Olive oil, native	17041404	1
Caelo	Olive oil, native	17051403	1
Caelo	Olive oil, native	18031506	1
Caelo	Olive oil, native	18041502	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Olive oil, native	23121413	1
Caelo	Olive oil, native	26071403	1
Caelo	Olive oil, native	26091403	1
Caelo	Olive oil, native	28051405	1
Bombastus	Olive oil, native	284143	2
Noweda	Olive oil, native	284143	1
S.Binder	Olive oil, native	284144	1
Fagron	Olive oil, native	284144	2
Bombastus	Olive oil, native	285160	1
Bombastus	Olive oil, native	285642	1
Gehe	Olive oil, native	289094	1
Gehe	Olive oil, native	289095	3
Noweda	Olive oil, native	295282	1
Caelo	Olive oil, native	295283	2
Caelo	Olive oil, native	297565	1
Anzag	Olive oil, native	297565	1
Bombastus	Olive oil, native	4613M-02000	1
Bombastus	Olive oil, native	4812A-02000	1
Sanacorp	Olive oil, native	5/AKD13	1
LAMOTTE/NOWEDA	Olive oil, native	5470202	1
Caelo	Olive oil, native	5470202	1
Lamotte/Phoenix	Olive oil, native	5623601	1
Gehe - Lamotte	Olive oil, native	5623601	1
Lamotte/Sanacorp	Olive oil, native	5470202	1
Caelo	Olive oil, native	5623601	1
Anzag/Lamotte	Olive oil, native	5623601	1
Gehe - Lamotte	Olive oil, native	5671601	2
Lamotte/Alliance	Olive oil, native	5623601	1
Caelo	Olive oil, native	5671601	1
Lamotte/Anzag	Olive oil, native	5671901	1
Henry Lamotte Oils GmbH...	Olive oil, native	5734701	1
Lamotte / Gehe	Olive oil, native	5734701	1
Lamotte/Anzag	Olive oil, native	5734701	1
Noweda	Olive oil, native	5813601	1
Henry Lamotte Oils GmbH	Olive oil, native	5813601	1
Lamotte/Alliance	Olive oil, native	5891202	1
Henry Lamotte/Anzag	Olive oil, native	590340003	1
Henry Lamotte Oils GmbH	Olive oil, native	5921100001	1
Henry Lamotte Oils GmbH	Olive oil, native	59211000001	1
noweda	Olive oil, native		1
Henry Lamotte Oils GmbH	Olive oil, native	8000652001	1
Lamotte/Kehr	Olive oil, native	8000821001	1
Henry Lamotte Oils	Olive oil, native	8001648001	1
Henry Lamotte Oils	Olive oil, native	8022366001	1
ACM	Olive oil, native	KNC2225J14	1
Caelo	Olive oil, native	CE834/ES	1
Bombastus	Olive oil, native	259542	1
Bombastus	Olive oil, native	265359	1
Sanacorp	Olive oil, native		1
Bombastus	Olive oil, native	256053	1
Gehe	Olive oil, native	265359	1
Bombastus	Olive oil, native	266513	1
Noweda	Olive oil, native	267376	1
Bombastus	Olive oil, native	267376	2
Bombastus	Olive oil, native	268118	1
Gehe	Olive oil, native	267376	1
Noweda	Olive oil, native	268118	1
Bombastus	Olive oil, native	269674	1
Gehe	Olive oil, native	269674	3
Bombastus	Olive oil, native	270971	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Olive oil, native	272270	2
Bombastus	Olive oil, native	272270	4
Phönix 09.11.12	Olive oil, native	272270	1
Sanacorp	Olive oil, native	273659	1
Bombastus	Olive oil, native	273659	3
Gehe	Olive oil, native	274236	1
Ichthyol/Phönix	Olive oil, native	274236	1
Caelo	Olive oil, native	1031302	1
Caelo	Olive oil, native	123773	1
Caelo	Olive oil, native	4041303	1
Caelo	Olive oil, native	4101312	1
Caelo	Olive oil, native	406212000	1
Caelo	Olive oil, native	503312000	1
Caelo	Olive oil, native	6111208	1
Caelo	Olive oil, native	7071206	1
Bombastus	Olive oil, native	8022014A	1
Caelo	Olive oil, native	713I-02000	1
Caelo	Olive oil, native	11063901	1
Bombastus	Olive oil, native	110816-1	1
Caelo	Olive oil, native	11063906	1
Caelo	Olive oil, native	10111201	1
Caelo	Olive oil, native	11227001	1
Caelo	Olive oil, native	11227002	2
Caelo	Olive oil, native	11227011	2
Kehr	Olive oil, native	11376906	1
Gehe	Olive oil, native	11C16-N08	1
Noweda	Olive oil, native	12077701	1
Anzag 31.10.12 ek: 2,78EUR	Olive oil, native	12077701	1
Cealo/Gehe	Olive oil, native	12077701	1
Caelo	Olive oil, native	120777027-13	1
Caelo	Olive oil, native	1207771	1
Resana GmbH	Olive oil, native	1208-3652	1
Euro OTC	Olive oil, native	1212022-01	2
Phönix 28.12.12	Olive oil, native	12141305	1
Caesar & Loretz GmbH, 4. . .	Olive oil, native	12141306	1
Cealo/Gehe	Olive oil, native	12141306	3
Sanacorp	Olive oil, native	12141306	2
Caelo	Olive oil, native	12141308	3
Noweda	Olive oil, native	12141309	1
Cae	Olive oil, native	12141309	1
Caelo	Olive oil, native	12141310	1
Caelo	Olive oil, native	1214306	1
Caelo	Olive oil, native	12274102	3
Phönix	Olive oil, native	12274106	1
Caelo	Olive oil, native	12274108	1
Wolf Handelsgesellschaft	Olive oil, native	128972ES	1
Sanacorp	Olive oil, native	12274108	1
Fagron	Olive oil, native	12A18-N08	2
Fagron	Olive oil, native	12I07-N07	3
Gehe	Olive oil, native	12I07-N07	1
Fagron	Olive oil, native	12I07-N07	1
Fagron	Olive oil, native	13010904	1
Hustadt Apotheke	Olive oil, native	13010904	1
Caelo	Olive oil, native	1301094	1
Caelo	Olive oil, native	13021410	1
Caelo	Olive oil, native	13081204	1
Alliance Healthcare	Olive oil, native	13224703	1
Fagron	Olive oil, native	13224704	2
Caelo	Olive oil, native	13224706	4
Caelo	Olive oil, native	1324704	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Olive oil, native	13306103	1
Caesar & Loretz GmbH	Olive oil, native	13350101	1
Phoenix	Olive oil, native	13350103	1
Fagron	Olive oil, native	13c08-N06	1
Caelo	Olive oil, native	13350106	4
J.H. Müller	Olive oil, native	140115	1
Phoenix	Olive oil, native	1406009-01	1
Caelo	Olive oil, native	14111306	1
Caelo	Olive oil, native	141182303	1
Caelo	Olive oil, native	1509212000	1
Caelo	Olive oil, native	20061207	1
Caelo	Olive oil, native	20270912-4	1
Caelo	Olive oil, native	2113Q02000	1
Caelo	Olive oil, native	24011408	1
Caelo	Olive oil, native	24011303	1
Caelo	Olive oil, native	265359	1
Gehe	Olive oil, native	275757	3
Alliance	Olive oil, native	277050	1
Gehe	Olive oil, native	278583	2
Phönix, 07.11.2011	Olive oil, native	278583	1
Gehe	Olive oil, native	281209	2
Gehe	Olive oil, native	281898	2
Gehe	Olive oil, native	284144	1
Bombastus	Olive oil, native	28963	1
Caelo	Olive oil, native	29041301	1
Caelo	Olive oil, native	3313I-02000	1
Caelo	Olive oil, native	3212I-02000	1
Caelo	Olive oil, native	5107101	2
Henry Lamotte Oils GmbH	Olive oil, native	5107101	1
Gehe	Olive oil, native	5107101	1
Henry Lamotte	Olive oil, native	5107101	1
Lamotte / Gehe	Olive oil, native	5127501	1
Henry Lamotte	Olive oil, native	5232201	1
HenryLamotte	Olive oil, native	5232201	2
Lamotte	Olive oil, native	5232201	3
Großhandel	Olive oil, native	5232201	2
H. Lamotte	Olive oil, native	5232201	1
H.Lamotte/Noweda	Olive oil, native	5232201	1
Henry, Noweda	Olive oil, native	5232201	1
Henry Lamotte / Phoenix	Olive oil, native	5232201	1
Lamotte/GH/11.09.13	Olive oil, native	5271504	1
Henry Lamotte	Olive oil, native	5271504	1
Henry Lamotte / Phoenix	Olive oil, native	5271504	1
HenryLamotte/Phoenix	Olive oil, native	5272401	4
Phönix/Lamotte	Olive oil, native	5272401	1
Lamotte Phönix	Olive oil, native	5272401	2
Henry Lamotte/Anzag	Olive oil, native	5272401	1
Großhandel	Olive oil, native	5272401	1
Lamotte / Sanacorp	Olive oil, native	5283901	1
Lamotte / Gehe	Olive oil, native	5283901	2
Henry Lamotte Oils GmbH	Olive oil, native	5372701	1
Lamotte / Gehe	Olive oil, native	5372701	2
Lamotte	Olive oil, native	5372701	1
AHD, 03.08.2011	Olive oil, native	5372701	1
Caelo	Olive oil, native	5464901	2
Lamotte/V.d.L.	Olive oil, native	5464901	1
Lamotte/Noweda	Olive oil, native	5464901	1
henry lamotte oils gmbh	Olive oil, native	5464901	1
Klenk	Olive oil, native	5470202	1
Großhandel	Olive oil, native	5470202	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Olive oil, native	60081503	1
Lamotte / Gehe	Olive oil, native	L901	1
Henry Lamotte Oils Gmbh...	Olive oil, native	5470202	1
Lamotte/Gehe	Omega-3-fish oil	124247	3
Henry Lamotte	Omega-3-fish oil	1000002451	1
Euro OTC	Omega-3-fish oil	1311028-01	1
Euro OTC	Omega-3-fish oil	1603021-01	1
Henry Lamotte Oils/Sana...	Omega-3-fish oil	123780	1
Lamotte/Sanacorp	Omega-3-fish oil	123883	1
Caelo	Almond oil, refined	465957	1
Lamotte über Sanacorp	Borage oil, refined, stabili...	125354	1
Henry Lamotte	Sesame oil, refined	124494	1
Caelo	Almond oil, refined	2101405	1
Caelo	Evening primrose oil, refined	3878-124466	1
Caelo	Almond oil, refined	6031506	1
Sanacorp	Wheatgerm oil, refined	1/AKD2	1
Lamotte/AHD	Evening primrose oil, refined	1000002386	1
Sanacorp	Sesame oil, refined	1000002705	1
Caelo	Almond oil, refined	1000003242	1
Henry Lamotte Oils / No...	Almond oil, refined	1000003242	1
Caelo	Evening primrose oil, refined	1000003443	1
Henry Lamotte Oils GmbH	Sesame oil, refined	1000003864	1
Henry Lamotte/Noweda	Almond oil, refined	1000004431	1
Henry Lamotte	Almond oil, refined	1000004431	1
Lamotte/Noweda	Almond oil, refined	1000004431	1
Bombastus	Corn oil, refined	10032015A	1
Caelo	Peanut oil, refined	1232304	1
Taoasis/Noweda	Almond oil, refined	130973	1
Caelo	Sesame oil, refined	13172908	2
Caelo	Sunflower oil, refined	13225102	2
Fagron	Almond oil, refined	13236411	1
Caelo	Wheatgerm oil, refined	13261603	7
Caelo	Wheatgerm oil, refined	13261604	3
Caelo	Apricot kernel oil, refined	13357604	1
Phönix	Peanut oil, refined	13358801	1
Caesar & Loretz GmbH	Peanut oil, refined	13358803	1
PHOENIX 15.05.2014	Borage oil, refined, stabili...	13372104	1
Fagron	Avocado oil, refined	13361406	1
Caelo	Sesame oil, refined	13397402	3
PHOENIX 08.05.2014	Evening primrose oil, refined	14002501	1
Caelo	Borage oil, refined, stabili...	13372110	1
Euro OTC	Evening primrose oil, refined	140133-01	1
Caelo	Wheatgerm oil, refined	14024501	6
Caelo	Wheatgerm oil, refined	14024502	1
PHOENIX 03.06.2014	Almond oil, refined	14042705	1
Caelo	Almond oil, refined	14042710	1
Caelo	Peanut oil, refined	14111401	1
Caesar & Loretz GmbH Hi...	Peanut oil, refined	14115404	1
Caesar & Loretz GmbH	Peanut oil, refined	14115404	1
Phoenix	Peanut oil, refined	14115405	1
Caelo	Almond oil, refined	141215804	1
Caelo	Grape seed oil, refined	14127902	4
Hedinger/ Noweda	Almond oil, refined	14139311	1
Caelo	Sunflower oil, refined	14147101	3
Caelo	Sesame oil, refined	14170001	1
Phoenix	Sesame oil, refined	14170004	1
14215814	Almond oil, refined	14215814	1
Schubert / Phoenix	Sunflower oil, refined	1423F0-2027	2
AHD	Peanut oil, refined	14246901	1
Caelo	Wheatgerm oil, refined	14272514	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sunflower oil, refined	14306701	1
Caelo	Sesame oil, refined	14351804	5
Sanacorp	Peanut oil, refined	14355904	1
Phoenix	Peanut oil, refined	14355905	1
Caesar&Loretz GmbH/Anzag	Peanut oil, refined	14355905	1
Euro OTC	Avocado oil, refined	150101501	1
Caelo	Almond oil, refined	15078401	2
Caelo	Almond oil, refined	15080602	1
Noweda	Evening primrose oil, refined	1509005-01	2
Caelo	Wheatgerm oil, refined	15111603	4
Caesar & Loretz/Noweda	Evening primrose oil, refined	1512005-01	1
Fischer/Gehe	Almond oil, refined	15127306	1
Caelo	Grape seed oil, refined	15144509	2
Caelo	Grape seed oil, refined	15144505	1
Caelo	Evening primrose oil, refined	15171108	1
Noweda	Almond oil, refined	15197404	1
Bombastus	Almond oil, refined	15197407	1
Bombastus	Almond oil, refined	15197409	1
Caesar & Loretz GmbH	Peanut oil, refined	15208506	1
Caelo	Soy oil, refined	15239310	1
Caelo	Sesame oil, refined	15239405	4
Caesar & Loretz/Noweda	Peanut oil, refined	15273402	1
Caesar/Noweda	Borage oil, refined, stabili...	15305003	1
Caelo	Sesame oil, refined	15374006	3
Caelo	Wheatgerm oil, refined	15374101	5
Caesar & Loretz/Noweda	Almond oil, refined	15409704	1
Caelo	Almond oil, refined	15409716	1
Caelo	Almond oil, refined	15806001	1
Caelo	Sunflower oil, refined	15422801	1
Caelo	Sesame oil, refined	15422706	1
apomix 004	Evening primrose oil, refined	1605027-01	1
Caelo	Evening primrose oil, refined	1605027-01	1
Caesar & Loretz/Noweda	Almond oil, refined	16071002	1
Caesar & Loretz/Noweda	Almond oil, refined	16071003	1
Phönix	Almond oil, refined	16071006	2
Caesar & Loretz	Almond oil, refined	16071006	1
Euro OTC	Almond oil, refined	16071006	1
Noweda	Peanut oil, refined	16071104	1
Caelo	Grape seed oil, refined	16092401	1
Caelo	Grape seed oil, refined	16092409	3
Bombastus	Almond oil, refined	16115303	1
Caelo	Wheatgerm oil, refined	16137404	5
Caesar & Loretz/Noweda	Peanut oil, refined	16249502	1
GEHE	Borage oil, refined, stabili...	16266011	1
Caelo	Sunflower oil, refined	16275801	1
Caelo	Corn oil, refined	16260309	1
Gehe	Evening primrose oil, refined	16281306	1
Caelo	Evening primrose oil, refined	16281315	1
Caelo	Evening primrose oil, refined	1628315	1
Caelo	Sesame oil, refined	16312507	1
Phoenix	Wheatgerm oil, refined	16321301	1
Caelo	Wheatgerm oil, refined	16321301	3
PHÖNIX	Almond oil, refined	16327302	1
Caesar & Loretz GmbH	Almond oil, refined	16327302	1
GEHE	Almond oil, refined	16327302	1
Phönix	Almond oil, refined	16327303	1
Gehe	Almond oil, refined	16327308	1
Phoenix	Almond oil, refined	16327304	1
Caesar&Loretz	Almond oil, refined	16327311	1
Phoenix	Almond oil, refined	163327303	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Grape seed oil, refined	16339703	2
Caelo	Grape seed oil, refined	16339706	2
Sanacorp	Peanut oil, refined	16348001	1
Phoenix	Peanut oil, refined	16348001	1
Euro OTC	Almond oil, refined	1701036-01	3
Phoenix	Evening primrose oil, refined	1703020-01	1
Caelo	Avocado oil, refined	17038804	1
Noweda	Evening primrose oil, refined	1704018	1
Caelo	Avocado oil, refined	17038811	1
Caelo	Evening primrose oil, refined	1704018-01	1
GEHE	Evening primrose oil, refined	1704018-01	1
GEHE	Almond oil, refined	17042203	1
Caelo	Almond oil, refined	17042206	1
GEHE	Peanut oil, refined	17042802	2
GEHE	Evening primrose oil, refined	1710011-01	1
Euro OTC	Evening primrose oil, refined	1710011-02	1
AHD	Peanut oil, refined	17169001	1
Aponorm/Wepa	Peanut oil, refined	17169018	1
GEHE	Almond oil, refined	17176501	1
Caesar & Loretz/Noweda	Almond oil, refined	17176502	1
Caelo	Grape seed oil, refined	17206507	2
Caelo	Almond oil, refined	17207905	2
Caelo	Sunflower oil, refined	17242104	1
Caelo	Sesame oil, refined	17332102	1
Caelo	Evening primrose oil, refined	18009106	3
AHD	Almond oil, refined	18015002	1
Caelo	Evening primrose oil, refined	18009110	1
GEHE	Almond oil, refined	18015007	1
Caesar & Loretz/ Noweda	Almond oil, refined	18015007	1
Caelo	Sunflower oil, refined	18027004	1
Eigenherstellung	Almond oil, refined	181114HML01	1
Eigenherstellung	Almond oil, refined	270215HJK01	1
Bombastus	Sunflower oil, refined	281796	2
Bombastus	Sunflower oil, refined	285712	2
Bombastus	Wheatgerm oil, refined	283449	1
Bombastus	Sunflower oil, refined	285243	1
Bombastus	Sunflower oil, refined	291278	1
Bombastus	Sunflower oil, refined	292956	1
Bombastus	Sunflower oil, refined	295818	2
Bombastus	Sunflower oil, refined	298487	1
Bombastus	Sunflower oil, refined	299712	1
Eigenherstellung	Almond oil, refined	300615HML01	1
Bombastus	Almond oil, refined	304024	1
Audor Pharma	Almond oil, refined	304023	1
MPK	Almond oil, refined	36480-127127	1
Henry Lamotte Oils/Jenn...	Sesame oil, refined	5514801	3
Lamotte	Almond oil, refined	5615001	2
lamotte noweda,	Almond oil, refined	5615001	1
Lamotte/Phoenix	Almond oil, refined	5615001	1
Lamotte / Noweda	Almond oil, refined	5615003	1
Phoenix/Lamotte	Sesame oil, refined	5634501	1
lamotte/ Phoenix	Sesame oil, refined	5744101	1
Henry Lamotte Oils GmbH	Peanut oil, refined	5699601	1
Henry Lamotte	Sunflower oil, refined	5738004	1
Lamotte/Sanacorp	Sesame oil, refined	5744101	1
Lamotte/Sanacorp	Almond oil, refined	5748401	1
Lamotte über Phoenix	Sesame oil, refined	5872001	1
Lamotte	Sesame oil, refined	5872901	1
Henry Lamotte/ NOweda	Sesame oil, refined	5872901	1
Noweda	Avocado oil, refined	713301202	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte/Phoenix	Almond oil, refined	5919600001	1
Henry Lamotte/Sanacorp	Sesame oil, refined	8000531001	1
Lamotte	Sesame oil, refined	8000531001	2
Henry Lamotte	Sesame oil, refined	8000531001	1
Lamotte über Noweda	Sesame oil, refined	8000682001	2
Lamotte über Phoenix	Almond oil, refined	8000725001	1
Henry Lamotte / Phoenix	Almond oil, refined	8000915002	1
Lamotte/Noweda	Almond oil, refined	8000915002	1
Lamotte	Sesame oil, refined	8001557002	2
Henry Lamotte/Kehr	Sesame oil, refined	8001557002	1
Henry Lamotte	Sesame oil, refined	8001557002	2
Lamotte	Almond oil, refined	80031670001	1
Henry Lamotte/Noweda	Almond oil, refined	8002675001	1
Lamotte	Almond oil, refined	8003167001	2
Fagron	Almond oil, refined	8003167001	1
Hennry Lamotte/Noweda	Almond oil, refined	8003503001	1
Lamotte	Almond oil, refined	8003885002	2
Lamotte/Phönix	Sesame oil, refined	84280	1
Anzag	Evening primrose oil, refined	AP141023N	1
Audor Pharma	Evening primrose oil, refined	AP562-124	1
Caelo	Evening primrose oil, refined	AP094-115N	1
Lamotte / Gehe	Avocado oil, refined	130193	1
Fagron	Avocado oil, refined	12F25-N12	1
Sanacorp	Avocado oil, refined	713301202	1
Henry Lamotte Oils	Borage oil, refined, stabili...	123118	1
Henry Lamotte	Borage oil, refined, stabili...	145678	2
Phönix	Borage oil, refined, stabili...	12161207	1
Caelo	Borage oil, refined, stabili...	13170603	1
Caelo	Borage oil, refined, stabili...	13372101	2
Caelo	Peanut oil, refined	8061303	1
Caelo	Peanut oil, refined	107063	1
Caelo	Peanut oil, refined	1224302	1
Caesar&Loretz GmbH/Anzag	Peanut oil, refined	12323204	1
Euro OTC	Peanut oil, refined	1301015	1
Caesar&Loretz GmbH/Anzag	Peanut oil, refined	13306106	1
Noweda	Peanut oil, refined	13358804	2
Caelo	Peanut oil, refined	13G30-N05	1
Caelo	Peanut oil, refined	28031304	1
Caelo	Peanut oil, refined	28021406	1
Caelo	Peanut oil, refined	16011403	1
Caelo	Peanut oil, refined	5314601	1
Bombastus	Almond oil, refined	273579	1
Caelo	Almond oil, refined	2021303	1
Eigenherstellung	Almond oil, refined	61113HML01	1
Caelo	Almond oil, refined	10131-01950	1
Caelo	Almond oil, refined	11358910	1
Caesar & Loretz/ Noweda	Almond oil, refined	12357101	2
GEHE	Almond oil, refined	13178003	1
Noweda	Almond oil, refined	13178011	2
Phönix 11.11.2013	Almond oil, refined	13236401	1
Anzag/Leyh Pharma am	Almond oil, refined	13236404	1
Caelo	Almond oil, refined	133500006	1
Caelo	Almond oil, refined	1335006	1
Caelo	Almond oil, refined	14139034	1
Eigenherstellung	Almond oil, refined	240513HML01	1
Bombastus	Almond oil, refined	274565	1
Phoenix/Lamotte	Almond oil, refined	5214206	1
Noweda/Lamotte Oils	Almond oil, refined	5214206	1
Lamotte/GEHE	Almond oil, refined	5269301	3
Lamotte/Noweda	Almond oil, refined	5438302	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix/Lamotte	Almond oil, refined	5438302	2
Lamotte	Almond oil, refined	548302	1
lamotte/phönix	Almond oil, refined	5438302	1
Taoasis	Almond oil, refined	5511701-113571	1
Caelo	Almond oil, refined	0-240513-1	1
Phönix	Evening primrose oil, refined	12325105	1
Fagron	Evening primrose oil, refined	13E06-No9	1
Caelo	Evening primrose oil, refined	1312816	1
Phoenix	Evening primrose oil, refined	14002503	1
Lamotte	Sesame oil, refined	123488	2
Caelo	Sesame oil, refined	11224108	3
Euro OTC	Sesame oil, refined	1205407	1
Caelo	Sesame oil, refined	12054501	1
Phönix 12.12.12	Sesame oil, refined	12054508	1
Caelo	Sesame oil, refined	12054508	4
Caelo	Sesame oil, refined	13011104	4
Caelo	Sesame oil, refined	13011107	8
Caelo	Sesame oil, refined	13011109	1
Caelo	Sesame oil, refined	130111107	1
Lamotte/noweda	Sesame oil, refined	170713a	1
Caelo	Sesame oil, refined	25051305	1
Henry Lamotte Oils	Sesame oil, refined	5319201	2
Henry Lamotte	Sesame oil, refined	5514801	1
Bombastus	Sunflower oil, refined	274107	2
Caelo	Sunflower oil, refined	12009804	3
Caelo	Grape seed oil, refined	12149204	4
Caelo	Grape seed oil, refined	12149205	1
Caelo	Walnut oil, refined	13338302	1
Caelo	Wheatgerm oil, refined	12220304	2
Caelo	Grape seed oil, refined	13065309	4
Caelo	Wheatgerm oil, refined	12220306	2
Caelo	Wheatgerm oil, refined	12220307	1
Caelo	Wheatgerm oil, refined	13011501	5
Caelo	Wheatgerm oil, refined	13011502	1
Caelo	Wheatgerm oil, refined	13011505	7
Caelo	Wheatgerm oil, refined	13011507	3
Euro OTC	Wheatgerm oil, refined	1305015-01	1
Caeöo	Wheatgerm oil, refined	13261601	1
Fagron	Wheatgerm oil, refined	13D22-N13	1
Caelo	Rapeseed oil, refined	15214403	1
Caelo	Wheatgerm oil, refined	W-240513-2	1
Caelo	Rapeseed oil, refined	15214407	1
Caelo	Rapeseed oil, refined	13011008	1
Caelo	Rapeseed oil, refined	15296106	3
Caelo	Rapeseed oil, refined	16185004	1
Caelo	Rapeseed oil, refined	17213901	2
Caelo	Rapeseed oil, refined	12054406	2
Fagron	Rapeseed oil, refined	13011003	1
Lamotte	Sesame oil, refined	123826	3
Lamotte	Sesame oil, refined	124494	2
Euro OTC	Sesame oil, refined	124494	1
Caelo	Sesame oil, refined	10011503	1
Caelo	Sesame oil, refined	11081513	1
Henry Lamotte Oils GmbH.	Sesame oil, refined	5639301	1
Henry Lamotte / Noweda	Sesame oil, refined	5639301	1
Lamotte über Phoenix	Sesame oil, refined	5639301	2
Lamotte	Sesame oil, refined	5744101	1
Lamotte	Sesame oil, refined	582001	1
Lamotte über Phoenix	Sesame oil, refined	5821001	1
Lamotte	Sesame oil, refined	5872001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte über Phoenix	Sesame oil, refined	5872901	2
Noweda/Henry Lamotte Oils	Sesame oil, refined	5872901	1
Henry Lamotte Oils GmbH	Sesame oil, refined	1000002705	1
Lamotte über Noweda	Sesame oil, refined	1000002705	1
Henry Lomotte/Noweda	Sesame oil, refined	1000002705	1
Henry Lamotte	Sesame oil, refined	5872901	1
Henry Lamotte Oils GmbH...	Sesame oil, refined	1000002705	1
Hedinger/Noweda	Sesame oil, refined	1000005858	1
Caelo	Sesame oil, refined	17177701	1
Caelo	Sesame oil, refined	17177703	1
Caesar und Loretz GmbH	Sesame oil, refined	16312505	1
Henry Lamotte	Sesame oil, refined	2705	1
Henry Lamotte Oils GmbH	Sesame oil, refined	122972	1
Gehe	Sesame oil, refined	122972	1
Gehe / 23.01.13 Henry L...	Sesame oil, refined	122972	1
Henry Lamotte/ 26.06.13...	Sesame oil, refined	123488	1
Henry Lamotte Oils GmbH	Sesame oil, refined	123488	1
Lamotte/ Noweda	Sesame oil, refined	123488	1
Henry Lamotte	Sesame oil, refined	123826	1
Henry Lamotte/Sanacorp	Sesame oil, refined	123826	1
Caelo	Sesame oil, refined	11224103	1
Henr	Sesame oil, refined	123826	1
Caelo	Sesame oil, refined	10241803	1
Caelo	Sesame oil, refined	123826	1
Caelo	Sesame oil, refined	1205408	1
Caelo	Sesame oil, refined	1301105	1
Euro OTC	Sesame oil, refined	13011107	1
Caelo	Sesame oil, refined	14101302	1
Caelo	Sesame oil, refined	17011405	1
Caelo	Sesame oil, refined	21051305	1
Caelo	Sesame oil, refined	30051201	1
Lamotte	Sesame oil, refined	4711A-02016	1
HenryLamotte/Phoenix	Sesame oil, refined	5135601	1
Henry Lamotte	Sesame oil, refined	5319201	1
Lamotte/Krieger	Sesame oil, refined	5319201	1
Lamotte	Sesame oil, refined	5514801	3
Euro OTC	Soy lecithin, liquid	1404022-01	2
Euro OTC	Soy lecithin, liquid	14404022-01	1
Euro OTC	Soy lecithin, liquid	1404022V-02	2
Euro OTC	Soy lecithin, liquid	1503038	1
Euro OTC	Soy lecithin, liquid	1204028-01	1
Gehe - Lamotte	Soy lecithin, liquid	1305019-01	1
Euro OTC	Soy lecithin, liquid	20061209	1
Euro OTC	Soy lecithin, liquid	1305019-01	1
Euro OTC	Soy lecithin, liquid	614E3336	1
Caelo	Soy oil, refined	14194107	1
Caelo	Soy oil, refined	16081309	1
Caelo	Sunflower oil, refined	16275804	1
Caelo	Sunflower oil, refined	17073005	1
Wolf Handelsgesellschaft	Sunflower oil, refined	1822313	1
Henry Lamotte/Sanacorp	Sunflower oil, refined	5738004	1
Henry Lamotte	Sunflower oil, refined	8001230001	1
Bombastus	Sunflower oil, refined	259035	1
Bombastus	Sunflower oil, refined	271457	1
Caelo	Sunflower oil, refined	11000605	2
Gehe	Sunflower oil, refined	11126606	1
J.H. Müller	Sunflower oil, refined	140229	1
Wolf Handelsgesellschaft	Sunflower oil, refined	1692305SK	1
Bombastus	Sunflower oil, refined	275830	1
Noweda / Henry Lamotte	Sunflower oil, refined	5422401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Lamotte/L3	Sunflower oil, refined	5422401	1
Caelo	Grape seed oil, refined	14214606	2
Caelo	Grape seed oil, refined	15144505	1
Caelo	Grape seed oil, refined	15294607	2
Fagron	Grape seed oil, refined	16092401	1
Caelo	Grape seed oil, refined	16092402	1
Caelo	Grape seed oil, refined	16092410	2
Caelo	Grape seed oil, refined	16339705	2
Caelo	Grape seed oil, refined	10222901	1
Caelo	Grape seed oil, refined	12149210	1
Caelo	Grape seed oil, refined	12149211	1
Caelo	Grape seed oil, refined	12149212	2
Noweda	Grape seed oil, refined	13065309	1
VDL;	Wheatgerm oil, refined	13261603	1
VDL;3,44EUR;22.05.15	Wheatgerm oil, refined	14356001	1
Gehe	Grape seed oil, refined	12149204	1
Caelo	Grape seed oil, refined	12149203	1
Holdermann	Grape seed oil, refined	12149212	1
Vögele	Walnut oil, refined	2017113718	1
Caelo	Walnut oil, refined	13202401	3
Caelo	Walnut oil, refined	12214902	1
Hepart AG	Walnut oil, refined	12214904	1
Caelo	Walnut oil, refined	12214904	1
Voegelé	Walnut oil, refined	2013113013	2
Caelo	Walnut oil, refined	24101309	1
Caelo	Wheatgerm oil, refined	5011502	1
Caelo	Wheatgerm oil, refined	13261601	5
Caelo	Wheatgerm oil, refined	14272511	2
Caelo	Wheatgerm oil, refined	14272513	2
Caelo	Wheatgerm oil, refined	14356001	5
Noweda	Wheatgerm oil, refined	14356003	1
Caelo	Wheatgerm oil, refined	14356003	1
Caelo	Wheatgerm oil, refined	14356006	3
Caelo	Wheatgerm oil, refined	14366001	1
Caelo	Wheatgerm oil, refined	15111602	3
Caelo	Wheatgerm oil, refined	15111605	1
Phoenix	Wheatgerm oil, refined	15374101	1
Caelo	Wheatgerm oil, refined	15374104	1
Caelo	Wheatgerm oil, refined	15406201	1
Caelo	Wheatgerm oil, refined	15406205	1
Caelo	Wheatgerm oil, refined	15406206	1
Phoenix	Wheatgerm oil, refined	16137403	1
phönix	Wheatgerm oil, refined	16321301	1
Phönix	Wheatgerm oil, refined	17216701	1
Bombastus	Wheatgerm oil, refined	286986	1
Caelo	Wheatgerm oil, refined	17290301	1
Bombastus	Wheatgerm oil, refined	289212	1
Caelo	Wheatgerm oil, refined	292513	1
Bombastus	Wheatgerm oil, refined	295374	1
Bombastus	Wheatgerm oil, refined	301383	1
phönix	Wheatgerm oil, refined		1
Caelo	Wheatgerm oil, refined	213E-02784	1
Primavera	Wheatgerm oil, refined	3606132	1
Caelo	Wheatgerm oil, refined	11227305	2
Ichthyl Phönix	Wheatgerm oil, refined	12220302	1
Großhandel	Wheatgerm oil, refined	12220302	1
Caelo	Wheatgerm oil, refined	12220305	1
Caelo	Wheatgerm oil, refined	13011503	2
Caelo	Wheatgerm oil, refined	13011508	2
Caelo	Wheatgerm oil, refined	13261605	1

- 83 058 spectra from 501 *Apo-Ident* customers from a total of 24 983 batches from a further 310 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Triglycerides* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Triglycerides* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	3340	0	32 204
Type B	0	4841	80	88 248
Type C	0	8214	80	82 766

The substance/substance group *Triglycerides* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9702 %)	100.0000 % (> 99.8204 %)
Type B	100.0000 % (> 99.9768 %)	98.3743 % (> 98.3134 %)
Type C	100.0000 % (> 99.6665 %)	99.0354 % (> 98.9993 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30647	31350	1.82	46.86
30679	31729	2.43	86.54
30754	31496	6.59	31.95
30840	33060	4.43	31.37
30855	31708	2.65	69.38
30863	31496	8.05	31.70
30873	31350	4.76	46.78

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30878	31361	5.52	69.94
30981	30981	0.00	65.40
31016	31016	0.00	55.36
31034	31523	0.94	31.27
31121	31121	0.00	61.60
31122	31122	0.00	20.87
31161	31161	0.00	61.65
31172	31172	0.00	75.12
31178	31729	4.87	84.32
31210	31210	0.00	57.44
31215	31215	0.00	47.86
31216	31216	0.00	64.85
31217	31217	0.00	53.35
31246	31246	0.00	32.52
31250	31250	0.00	29.14
31274	31274	0.00	52.16
31322	31322	0.00	42.47
31326	33064	4.96	80.09
31332	31332	0.00	49.75
31336	31336	0.00	34.60
31337	31770	2.05	36.73
31344	31274	2.22	52.57
31346	31801	1.79	33.07
31350	31350	0.00	45.28
31356	32195	1.41	34.32
31361	31361	0.00	65.78
31362	30981	4.63	62.64
31374	31274	2.11	51.96
31375	31375	0.00	28.03
31376	31350	2.82	45.12
31378	31491	4.23	55.27
31380	31216	2.23	64.14
31382	31382	0.00	41.38
31384	31161	2.94	61.90
31460	31460	0.00	66.30
31491	31491	0.00	54.03
31496	31496	0.00	31.78
31523	31523	0.00	28.40
31540	31729	2.93	86.37
31572	31572	0.00	65.31
31584	31382	1.22	42.13
31585	31801	3.33	35.70
31589	31572	2.79	66.38
31609	31609	0.00	51.00
31662	31801	2.39	32.29
31677	31361	6.92	67.90
31708	31708	0.00	65.86
31729	31729	0.00	88.27
31770	31770	0.00	36.54
31775	31775	0.00	67.60
31792	31729	0.96	87.57
31801	31801	0.00	33.80
31825	31361	5.49	64.75

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31833	31833	0.00	16.51
31834	31016	4.72	57.85
31875	31729	7.75	82.94
32195	32195	0.00	32.81
32680	32680	0.00	33.41
32683	32683	0.00	28.78
32730	31210	2.16	56.22
32745	31217	3.20	52.31
32989	30981	6.30	60.29
32991	31801	3.82	34.23
33059	33059	0.00	30.50
33060	33060	0.00	31.91
33064	33064	0.00	73.53
33067	33067	0.00	32.00
33070	33070	0.00	38.07
33254	31322	6.30	41.56
33255	31122	2.67	19.07
33309	33309	0.00	59.67
33313	33313	0.00	46.12
33510	30981	4.51	63.00
33876	31322	4.34	40.66
33916	33916	0.00	65.67
33917	33917	0.00	50.13

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Triglycerides, medium-chain**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 31483-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Triglycerides, medium-chain; Oleum neutrale; Triglycerida saturada media

Special notes

When selecting the *Triglycerides, medium-chain* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Triglycerides, medium-chain*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Triglycerides, m...	13392404	31483	40	AR-14-FG-012976-01
Caelo	Triglycerides, m...	13348507	31484	40	20131028
Euro OTC	Triglycerides, m...	1402019-01	31485	40	20140317
Euro OTC	Triglycerides, m...	1509022	31771	40	20151007*
Euro OTC	Triglycerides, m...	1511011	31804	40	20151125*
Euro OTC	Triglycerides, m...	1601050	31955	40	20160218*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 240 spectra of 6 reference samples from the substance/substance group *Triglycerides, medium-chain*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 6 different batches.
- 35 304 spectra from a total of 834 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 560 spectra from 17 reference samples of the substance/substance group *Triglycerides, medium-chain*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
zu erfassen	Triglycerides, medium-chain	zu erfassen	30174	40
Caelo	Triglycerides, medium-chain	13392404	31483 [†]	20
Caelo	Triglycerides, medium-chain	13348507	31484 [†]	20
Euro OTC	Triglycerides, medium-chain	1402019-01	31485 [†]	20
Euro OTC	Triglycerides, medium-chain	1509022	31771 [†]	20
Euro OTC	Triglycerides, medium-chain	1511011	31804 [†]	20
Euro OTC	Triglycerides, medium-chain	1601050	31955 [†]	20
Caelo	Triglycerides, medium-chain	161243	32204	40
Euro OTC	Triglycerides, medium-chain	1603041	32278	40
Euro OTC	Triglycerides, medium-chain	1606035	32464	40
Caelo	Triglycerides, medium-chain	162385	32697	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Triglycerides, medium-chain	1701023	32983	40
Caelo	Triglycerides, medium-chain	170160	33077	40
Caelo	Triglycerides, medium-chain	170784	33312	40
Caelo	Triglycerides, medium-chain	172177	33672	40
Euro OTC	Triglycerides, medium-chain	1709014	33845	40
Euro OTC	Triglycerides, medium-chain	1801034	33997	40

- 93314 spectra from a total of 2151 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1242 spectra from 352 *Apo-Ident* customers from 314 batches from the substance/substance group *Triglycerides, medium-chain*.
- These include spectra of independent samples from 302 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Triglycerides, medium-chain	109948	1
THC Pharma	Triglycerides, medium-chain	116690	1
Bionorica ethics	Triglycerides, medium-chain	170517	1
Bionorica/Phoenix	Triglycerides, medium-chain	171102	1
Caelo	Triglycerides, medium-chain	17194902	1
Caelo	Triglycerides, medium-chain	17236202	6
Caelo	Triglycerides, medium-chain	17361001	10
Bombastus	Triglycerides, medium-chain	17361001	1
Bionorica Ethics/Noweda	Triglycerides, medium-chain	85113	3
Bionorica ethics	Triglycerides, medium-chain	91713	2
Bionorica	Triglycerides, medium-chain	91719	2
Sanacorp	Triglycerides, medium-chain	93574	1
Bionorica/Krieger	Triglycerides, medium-chain	98569	1
THC/Noweda	Triglycerides, medium-chain	93572	1
bionorica	Triglycerides, medium-chain	93574	1
GEHE	Triglycerides, medium-chain	98569	1
Bionorica ethics/Noweda	Triglycerides, medium-chain	105039	2
THC	Triglycerides, medium-chain	105041	2
Sanacorp WE:27.09.17 ...	Triglycerides, medium-chain	107077	1
Bionorica Ethics	Triglycerides, medium-chain	109709	1
Bionorica/Kehr	Triglycerides, medium-chain	109951	1
THC	Triglycerides, medium-chain	105043	1
Caelo	Triglycerides, medium-chain	11301106	14
Caelo	Triglycerides, medium-chain	13159607	1
Caelo	Triglycerides, medium-chain	13348501	5

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triglycerides, medium-chain	14098306	4
Caelo	Triglycerides, medium-chain	14219003	11
Sanacorp	Triglycerides, medium-chain	14219003	3
Caelo	Triglycerides, medium-chain	14219007	10
Noweda	Triglycerides, medium-chain	14219007	1
Caelo	Triglycerides, medium-chain	14219009	8
Caelo	Triglycerides, medium-chain	14380802	9
Caelo	Triglycerides, medium-chain	14380803	9
Caelo	Triglycerides, medium-chain	14380804	20
Fagron	Triglycerides, medium-chain	14G25B03	1
Fagron	Triglycerides, medium-chain	14G25-B03	3
Fagron	Triglycerides, medium-chain	14G25-B03-296937	5
Fagron	Triglycerides, medium-chain	14G25-B03-302294	10
Euro OTC	Triglycerides, medium-chain	1504029-02	9
Euro OTC	Triglycerides, medium-chain	1509022-01	4
Bionorica ethics/Sanaco...	Triglycerides, medium-chain	150930	1
Caelo	Triglycerides, medium-chain	15139604	6
Caelo	Triglycerides, medium-chain	15139605	11
Phoenix	Triglycerides, medium-chain	15139605	1
Miglyol	Triglycerides, medium-chain	15139605	1
Caelo	Triglycerides, medium-chain	15139607	17
Caelo	Triglycerides, medium-chain	15139608	25
Noweda	Triglycerides, medium-chain	15139608	1
Sanacorp	Triglycerides, medium-chain	15139608	1
Caelo	Triglycerides, medium-chain	15270004	1
Caelo	Triglycerides, medium-chain	15272003	26
Caelo	Triglycerides, medium-chain	15272004	32
/Sanacorp 02.05.2016	Triglycerides, medium-chain	15272004	1
Caelo	Triglycerides, medium-chain	15272005	11
Caelo	Triglycerides, medium-chain	15272006	13
Caelo	Triglycerides, medium-chain	15272007	6
Caelo	Triglycerides, medium-chain	152722007	1
Caelo	Triglycerides, medium-chain		4
Caelo	Triglycerides, medium-chain	15272008	2
Caelo	Triglycerides, medium-chain	15379602/170314M	1
Henry Lamotte Oils GmbH...	Triglycerides, medium-chain	15379602	1
Caelo	Triglycerides, medium-chain	15379608	17
Sanacorp	Triglycerides, medium-chain	15379608	1
Caolo/Noweda	Triglycerides, medium-chain	15379608	1
Caelo	Triglycerides, medium-chain	15379609	9
Caelo	Triglycerides, medium-chain	15379610	9
Caesar & Loretz GmbH Hi...	Triglycerides, medium-chain	15379610	1
Euro OTC	Triglycerides, medium-chain	15379610	1
Hubertus Apotheke	Triglycerides, medium-chain	15379610	1
Caelo	Triglycerides, medium-chain	15379611	2
Ebert 12.10.16	Triglycerides, medium-chain	15379611	1
Caelo	Triglycerides, medium-chain	15379612	5
Euro OTC	Triglycerides, medium-chain	15379612	1
Caelo	Triglycerides, medium-chain	1572006	1
Fagron	Triglycerides, medium-chain	15H03-B03	2
Fagron	Triglycerides, medium-chain	15H03-B03-318529	5
Fagron	Triglycerides, medium-chain	15H03-B03318529	1
Fagron	Triglycerides, medium-chain	15H03-B03-318530	1
Fagron	Triglycerides, medium-chain	15H03-B03-311235	1
Euro OTC	Triglycerides, medium-chain	1606035-01	13
Fagron	Triglycerides, medium-chain	15H03-B03-318529	1
Caelo	Triglycerides, medium-chain	161243047	1
Caelo	Triglycerides, medium-chain	16124304	3
Euro OTC	Triglycerides, medium-chain	1601050-01	27
Caelo	Triglycerides, medium-chain	16124306/170518M	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triglycerides, medium-chain	16124307	30
Caelo	Triglycerides, medium-chain	16124308	8
Caelo	Triglycerides, medium-chain	16124309	1
Caelo	Triglycerides, medium-chain	17016001	12
Caesar & Loretz GmbH, P...	Triglycerides, medium-chain	17016001	1
Caelo	Triglycerides, medium-chain	17016002	12
Caesar und Loretz GmbH	Triglycerides, medium-chain	17016007	1
Caelo	Triglycerides, medium-chain	17016007	6
Caelo	Triglycerides, medium-chain	17194901	2
Caelo	Triglycerides, medium-chain	17016008	1
Caelo	Triglycerides, medium-chain	170420MT/16124308	1
Caelo	Triglycerides, medium-chain	171123M/17016009	1
Caelo	Triglycerides, medium-chain	1706007	1
Caelo	Triglycerides, medium-chain	180109M/17236204	1
Caelo	Triglycerides, medium-chain	180206M/17236204	1
THC Pharm	Triglycerides, medium-chain		1
THC Pharma/Phönix	Triglycerides, medium-chain		1
Caelo	Triglycerides, medium-chain	2014030260	1
THC/Sana	Triglycerides, medium-chain	2015070139	2
Sanacorp	Triglycerides, medium-chain	6/AKD12	1
Sanacorp	Triglycerides, medium-chain	58/AKD4	1
THC Pharma/Phönix	Triglycerides, medium-chain	90059	2
Caelo	Triglycerides, medium-chain	6061201	1
v.d.l	Triglycerides, medium-chain	10356201	1
THC/Noweda	Triglycerides, medium-chain	91714	1
Caelo	Triglycerides, medium-chain	10356205	1
Fagron	Triglycerides, medium-chain	110321	1
Caelo	Triglycerides, medium-chain	11207502	2
Noweda	Triglycerides, medium-chain	11207501	1
Caelo	Triglycerides, medium-chain	11301105	1
Caelo	Triglycerides, medium-chain	11207501	1
Fiebig	Triglycerides, medium-chain	11207502	1
Phönix	Triglycerides, medium-chain	11301106	1
Fagron	Triglycerides, medium-chain	11L19-N20	1
Caelo	Triglycerides, medium-chain	12078603	14
Caelo	Triglycerides, medium-chain	13081205	1
Fagron	Triglycerides, medium-chain	12A20-N16	2
Caelo	Triglycerides, medium-chain	27101202	1
Caelo	Triglycerides, medium-chain	6051503	1
Bionorica	Triglycerides, medium-chain	87113	1
Caelo	Triglycerides, medium-chain	11031502	1
AHD 23.12.2014	Triglycerides, medium-chain	14219003	1
Gehe	Triglycerides, medium-chain	14D23-B02-293330	1
Sanacorp	Triglycerides, medium-chain	56/AKD3	1
Sanacorp/Bionorica Ethics	Triglycerides, medium-chain	83204	1
Phönix 07.11.12	Triglycerides, medium-chain	58447	1
Bionorica	Triglycerides, medium-chain	58448	11
Bionorica Ethics	Triglycerides, medium-chain	62424	1
Bionorica	Triglycerides, medium-chain	79406	1
Caelo	Triglycerides, medium-chain	8011301	1
Caelo	Triglycerides, medium-chain	1115805	1
Caelo	Triglycerides, medium-chain	2101305	1
Fagron	Triglycerides, medium-chain	101001	1
Apotheke in Steinkuhl	Triglycerides, medium-chain	101110	1
Caelo	Triglycerides, medium-chain	11154601	1
Caelo	Triglycerides, medium-chain	11207503	1
Caelo	Triglycerides, medium-chain	11301109	3
Caesar & Loretz GmbH	Triglycerides, medium-chain	12078603	1
Caelo	Triglycerides, medium-chain	12078604	12
Sanacorp	Triglycerides, medium-chain	12078604	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
gehe, 23.01.2013	Triglycerides, medium-chain	12078604	1
Caelo	Triglycerides, medium-chain	12078605	4
Euro OTC	Triglycerides, medium-chain	1210002-01	5
Caelo	Triglycerides, medium-chain	12101303	1
Caelo	Triglycerides, medium-chain	1225604	1
Caelo	Triglycerides, medium-chain	1225605	1
Caelo	Triglycerides, medium-chain	12256504	6
Caelo	Triglycerides, medium-chain	12256505	25
Noveda	Triglycerides, medium-chain	12256505	1
Phönix	Triglycerides, medium-chain	12256505	1
Gehe	Triglycerides, medium-chain	12256505	1
Caelo	Triglycerides, medium-chain	1226505	1
Gehe 06.05.2013	Triglycerides, medium-chain	12256505	1
Caelo	Triglycerides, medium-chain	1235601	1
Caelo	Triglycerides, medium-chain	12356901	6
Caelo	Triglycerides, medium-chain	12356903	11
Phönix	Triglycerides, medium-chain	12356901	2
Caelo	Triglycerides, medium-chain	12356904	1
Caelo	Triglycerides, medium-chain	12356908	26
Caelo	Triglycerides, medium-chain	12356909	4
Caelo	Triglycerides, medium-chain	1256505	1
Fagron	Triglycerides, medium-chain	12I04-N07	9
Fagron	Triglycerides, medium-chain	12K07-N22	5
Caelo	Triglycerides, medium-chain	13095304	11
Caelo	Triglycerides, medium-chain	13095306	3
celo/phönix	Triglycerides, medium-chain	13095306	1
Caelo	Triglycerides, medium-chain	13098601	16
Großhandel	Triglycerides, medium-chain	13098601	1
Ebert	Triglycerides, medium-chain	13098601	1
Caelo	Triglycerides, medium-chain	13098605	7
Caelo	Triglycerides, medium-chain	13177101	9
Gehe	Triglycerides, medium-chain	13177102	3
Phönix	Triglycerides, medium-chain	13177102	1
Caelo	Triglycerides, medium-chain	13177102	24
Caesar&lorenz,Phoenix	Triglycerides, medium-chain	13177102	1
Caelo	Triglycerides, medium-chain	13177111	1
Caelo	Triglycerides, medium-chain	13348502	7
Sanacorp	Triglycerides, medium-chain	13348502	2
Caelo	Triglycerides, medium-chain	13348503	8
Caelo	Triglycerides, medium-chain	13348506	5
Caelo	Triglycerides, medium-chain	13348507	1
Fagron	Triglycerides, medium-chain	13E24-N04	14
Fagron	Triglycerides, medium-chain	13H02N03	1
Caelo	Triglycerides, medium-chain	14098305	1
Caelo	Triglycerides, medium-chain	1803111	2
Caelo	Triglycerides, medium-chain	2113Q01998	1
Caelo	Triglycerides, medium-chain	20270312-3	1
Caelo	Triglycerides, medium-chain	18111309	1
HCT Pharma	Triglycerides, medium-chain	2012100833	1
Caelo	Triglycerides, medium-chain	4110i-01998	1
Caelo	Triglycerides, medium-chain	31051303	1
Caelo	Triglycerides, medium-chain	18031307	1
Euro OTC	Triglycerides, medium-chain	58448	2
Euro OTC	Triglycerides, medium-chain	67940	2
Bionorica	Triglycerides, medium-chain	69719	2
Bionorica Ethics	Triglycerides, medium-chain	75316	1
Bionorica/Sanacorp WE:	Triglycerides, medium-chain	85532	1
EuRho	Triglycerides, medium-chain	4121201	1
Euro OTC	Triglycerides, medium-chain	1105005-01	2
Euro OTC	Triglycerides, medium-chain	1111037-01	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
EuRho	Triglycerides, medium-chain	1111037-01	1
Caelo	Triglycerides, medium-chain	614M1998	1
Phönix Gotha	Triglycerides, medium-chain	1111037-01	1
Caesar/Noweda	Triglycerides, medium-chain	11301104	1
anzag, 20.11.2012	Triglycerides, medium-chain	11L19-N20	2
Phönix Gotha	Triglycerides, medium-chain	1205015-01	2
Euro OTC	Triglycerides, medium-chain	1205015-01	3
Fagron	Triglycerides, medium-chain	1210002-01	3
Phönix Gotha	Triglycerides, medium-chain	1210002-01	2
Phönix	Triglycerides, medium-chain	1210002-01	2
Caelo	Triglycerides, medium-chain	1210002-01	4
phönix 22.01.2013	Triglycerides, medium-chain	1210002-01	1
Phönix Gotha	Triglycerides, medium-chain	1302004-01	3
Euro OTC	Triglycerides, medium-chain	130200401	1
Euro OTC	Triglycerides, medium-chain	1302004-01	13
Euro OTC	Triglycerides, medium-chain	1302004-02	8
Caelo	Triglycerides, medium-chain	13095303	2
Caelo	Triglycerides, medium-chain	13124310	1
Caesar/Noweda	Triglycerides, medium-chain	13177109	1
Caelo	Triglycerides, medium-chain	13392404	2
Hansen&Rosenthal/Phönix	Triglycerides, medium-chain	1406016-01	1
Euro OTC	Triglycerides, medium-chain	1406016-01	1
Phönix Gotha	Triglycerides, medium-chain	1411009-01	2
Euro OTC	Triglycerides, medium-chain	1411009-02	6
Phönix Gotha	Triglycerides, medium-chain	1411009-02	1
Caelo	Triglycerides, medium-chain	1411009-02	1
Caelo	Triglycerides, medium-chain	14330503	1
Caelo	Triglycerides, medium-chain	15012016	1
Phönix Gotha	Triglycerides, medium-chain	1504029-01	2
Caelo	Triglycerides, medium-chain	150402902	1
Euro OTC	Triglycerides, medium-chain	1509022-01G	1
Euro OTC	Triglycerides, medium-chain	1511011-01	9
Caelo	Triglycerides, medium-chain	15239802	1
Caelo	Triglycerides, medium-chain	1537908	1
Caelo	Triglycerides, medium-chain	15H03-B03-312351	1
Fagron	Triglycerides, medium-chain	15H03-B03-328583	2
Gehe	Triglycerides, medium-chain	1601050-0	1
Caelo	Triglycerides, medium-chain	1603041-02	3
Euro OTC	Triglycerides, medium-chain	1603041-02	14
Ichthyol/Phönix	Triglycerides, medium-chain	1603041-02	2
Euro OTC	Triglycerides, medium-chain	1606035-02	11
Euro OTC	Triglycerides, medium-chain	1606035-03	4
Euro OTC	Triglycerides, medium-chain	1701023-01	1
Euro OTC	Triglycerides, medium-chain	1709014-01	3
Caelo	Triglycerides, medium-chain	18118804	1
THC/Phönix	Triglycerides, medium-chain	2014030260	2
Caelo	Triglycerides, medium-chain	50775165	1
Caelo	Triglycerides, medium-chain	3413I-01998	1
Pharimex	Triglycerides, medium-chain	75316	2
THC Pharma/Phönix	Triglycerides, medium-chain	85533	1
THC Pharma/Phönix	Triglycerides, medium-chain	93571	1
ACM	Triglycerides, medium-chain	KNC0254J15	5
Euro OTC	Triglycerides, medium-chain	1411009-01	6
Caelo	Triglycerides, medium-chain	1539605	1
Fagron	Triglycerides, medium-chain	15H03-B03-312351	9
Fagron	Triglycerides, medium-chain	15H03-B03-525120	8
Euro OTC	Triglycerides, medium-chain	2614A-01998	1
Sanacorp	Triglycerides, medium-chain	12356908	1
Bionorica Ethics/AHD	Triglycerides, medium-chain	62424	2
Sanacorp	Triglycerides, medium-chain	12356903	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
BionoricPhoenix	Triglycerides, medium-chain	958448	1
THC Pharm/Sanacorp	Triglycerides, medium-chain	91714	1
Bionorica	Triglycerides, medium-chain	87112	2
Bionorica/Noweda	Triglycerides, medium-chain	93566	1
Bionorica ethics/Sanaco...	Triglycerides, medium-chain	93569	1
THC Pharm	Triglycerides, medium-chain	93571	1
THC Pharm / Gehe	Triglycerides, medium-chain	93573	2
THC Pharma/Noweda	Triglycerides, medium-chain	93573	1
Bionorica Ethics	Triglycerides, medium-chain	94875	1
GEHE	Triglycerides, medium-chain	95080	1
Bionorica ethics	Triglycerides, medium-chain	94875	1
Bionorica	Triglycerides, medium-chain	95080	1
THC Pharm / Phönix	Triglycerides, medium-chain	95083	1
Bionorica Ethics/Noweda	Triglycerides, medium-chain	95112	1
Bionorica ethics/Sanaco...	Triglycerides, medium-chain	97884	1
Bionorica/Noweda	Triglycerides, medium-chain	98573	1
Bionorica/Noweda	Triglycerides, medium-chain	103899	2
Bionorica/Noweda	Triglycerides, medium-chain	103900	1
Bionorica/Noweda	Triglycerides, medium-chain	105038	1
bio	Triglycerides, medium-chain	105038	1
Bionorica Ethics	Triglycerides, medium-chain	105039	1
GEHE	Triglycerides, medium-chain	105039	1
GEHE	Triglycerides, medium-chain	107077	2
GEHE	Triglycerides, medium-chain	107078	1
Bionorica/Noweda	Triglycerides, medium-chain	107080	1
Bionorica Dronabinol-Set	Triglycerides, medium-chain	107078	1
Bionorica/ Noweda	Triglycerides, medium-chain	109950	1
GEHE	Triglycerides, medium-chain	109949	1
Gehe	Triglycerides, medium-chain	109951	1
Fagron	Triglycerides, medium-chain	116638	1
Bionorica/Kehr	Triglycerides, medium-chain	116686	1
GEHE	Triglycerides, medium-chain	116687	1
Bionorica Ethics	Triglycerides, medium-chain	109951	1
Fagron	Triglycerides, medium-chain	12/04-N07	1
Fagron	Triglycerides, medium-chain	120604	2
bionoricaEthiks/Phoenix	Triglycerides, medium-chain	116687	1
Fagron	Triglycerides, medium-chain	140514	1
EuRho	Triglycerides, medium-chain	1406016-01	1
	Triglycerides, medium-chain	1406016-01	1
Euro OTC	Triglycerides, medium-chain	1411009-002	1
EuRho/Phönix	Triglycerides, medium-chain	1411009-02	1
Euro OTC	Triglycerides, medium-chain	141100902	1
Caesar/Alliance	Triglycerides, medium-chain	14219007	1
Fagron	Triglycerides, medium-chain	14D08-N01	2
Fagron	Triglycerides, medium-chain	14D11-N01	1
Fagron	Triglycerides, medium-chain	14G25-B03-305014	1
Fagron	Triglycerides, medium-chain	14K28-N01	2
Euro OTC	Triglycerides, medium-chain	1504029-01	7
Fagron	Triglycerides, medium-chain		2
Euro OTC	Triglycerides, medium-chain	150402901	1
Anzag/Henry Lamotte	Triglycerides, medium-chain	1504029-02	1
Caelo	Triglycerides, medium-chain	1504029-02	1
Fargon/Phönix	Triglycerides, medium-chain	150518	1
Fagron	Triglycerides, medium-chain	150518	1
Caelo	Triglycerides, medium-chain	150518	1
Caelo	Triglycerides, medium-chain	15080603	1
Hansen & Rosenthal/Noweda	Triglycerides, medium-chain	1509022-01	1
Euro OTC	Triglycerides, medium-chain	151101101	1
Bionorica/Kehr	Triglycerides, medium-chain	151130	1
Bombastus	Triglycerides, medium-chain	15139608	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triglycerides, medium-chain	1527003	1
WE:05.04.2016 Sanacorp	Triglycerides, medium-chain	15272006	1
Caelo	Triglycerides, medium-chain	1537961	1
Fagron	Triglycerides, medium-chain	15403-B03-318529	1
Fagron	Triglycerides, medium-chain	15D24-B03	1
Fagron	Triglycerides, medium-chain	15D24-B03-307759	1
Caelo	Triglycerides, medium-chain	15H03-B03(BTM-Set)	1
Fagron	Triglycerides, medium-chain	15D24-B03-525046	1
Fagron	Triglycerides, medium-chain	15H03-B03-312924	5
Fagron	Triglycerides, medium-chain	15H03-B03-325574	2
Cremer Oleo GmbH & Co K...	Triglycerides, medium-chain	15H03-B03-525120	1
Fagron	Triglycerides, medium-chain	15H03-B03-525189	1
Fagron	Triglycerides, medium-chain	15H03-B03-525271	2
Phoenix	Triglycerides, medium-chain	15H03-B03-525271	1
Fagron	Triglycerides, medium-chain	15H03-B03-525225	1
Fagron	Triglycerides, medium-chain	15H03-B03-525383	1
Caelo	Triglycerides, medium-chain	15H-B03-525331	1
Caelo	Triglycerides, medium-chain	15H-B03-525120	1
Euro OTC	Triglycerides, medium-chain	160105001	2
Caelo	Triglycerides, medium-chain	1601050-01	2
apo-ident/ Phönix	Triglycerides, medium-chain	1601050-01	1
Gehe	Triglycerides, medium-chain	1601050-01	1
Euro OTC	Triglycerides, medium-chain	1603041-01	2
Euro OTC	Triglycerides, medium-chain	16030410	1
Gehe	Triglycerides, medium-chain	1603041-02	1
GEHE	Triglycerides, medium-chain	1603041-01	1
GEHE	Triglycerides, medium-chain	1603041-02	1
phönix	Triglycerides, medium-chain	1606035	1
hochhaus-Apotheke	Triglycerides, medium-chain	1606035-01	1
Pharmachem	Triglycerides, medium-chain	1606035-02	1
Caelo	Triglycerides, medium-chain	1606035-02	1
Apo-Ident/Noweda	Triglycerides, medium-chain	1606035-02	1
GEHE	Triglycerides, medium-chain	1606035-03	1
Gehe	Triglycerides, medium-chain	16124307	2
Caelo	Triglycerides, medium-chain	16124310	1
Gehe	Triglycerides, medium-chain	16124310	1
Caelo	Triglycerides, medium-chain	16288504	1
Fagron	Triglycerides, medium-chain	16H18-B01	2
Fagron	Triglycerides, medium-chain	16H18-B01-525420	3
Fagron	Triglycerides, medium-chain	16H18-B01-525477	1
Gehe	Triglycerides, medium-chain	16H18-B01-525477	1
Caelo	Triglycerides, medium-chain	16H18-B01-525506	1
Fagron	Triglycerides, medium-chain	16H18-B01-525531	1
GEHE	Triglycerides, medium-chain	17016001	1
Noweda	Triglycerides, medium-chain	17016001	1
Fiebig	Triglycerides, medium-chain	17016007	1
Gehe	Triglycerides, medium-chain	1709014-01	1
Defektur	Triglycerides, medium-chain	1710910-b	1
eigene Herstellung	Triglycerides, medium-chain	1712907-a	1
Noweda	Triglycerides, medium-chain	17236202	1
Fagron	Triglycerides, medium-chain	17F16-B02-341874	1
Fagron	Triglycerides, medium-chain	17F16-B02	1
Fagron	Triglycerides, medium-chain	17F16-B02-525622	1
GEHE	Triglycerides, medium-chain	1801034-01	1
THC Pharma/Noweda	Triglycerides, medium-chain	2014030260	1
Caelo	Triglycerides, medium-chain	59/AKD7	1
Euro OTC	Triglycerides, medium-chain	41141312	1
THC Pharma/Noweda	Triglycerides, medium-chain	2015070139	1
THC	Triglycerides, medium-chain	180522PMT0000105043	1
Euro OTC	Triglycerides, medium-chain	700282	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Triglycerides, medium-chain	61140578	1
Sanacorp	Triglycerides, medium-chain		1
Anzag	Triglycerides, medium-chain	12078604	2
Anzag	Triglycerides, medium-chain	12256505	1
Gehe	Triglycerides, medium-chain	1302004-01	2
Caelo	Triglycerides, medium-chain	613A-01998	1
Euro OTC	Triglycerides, medium-chain	1007009-01	1
Caelo	Triglycerides, medium-chain	101001	1
Caelo	Triglycerides, medium-chain	10356208	1
Noweda	Triglycerides, medium-chain	101001	1
Fagron	Triglycerides, medium-chain	10K09-0031	1
Caelo	Triglycerides, medium-chain	11154602	2
Anzag 11.10.12	Triglycerides, medium-chain	11301109	1
Phönix	Triglycerides, medium-chain	11320508	1
Fagron	Triglycerides, medium-chain	11C22-N03	1
Caelo	Triglycerides, medium-chain	11D29-N01	1
Fagron	Triglycerides, medium-chain	11G28-N07	1
Eurho	Triglycerides, medium-chain	1205015-01	1
Gehe	Triglycerides, medium-chain	12078603	1
Phönix 29.10.12	Triglycerides, medium-chain	12078604	1
Caelo	Triglycerides, medium-chain	12253802	5
Fagron	Triglycerides, medium-chain	12C27-N02	1
Fagron	Triglycerides, medium-chain	12I04-N07	1
Phagron GmbH / Anzag	Triglycerides, medium-chain	12K07-N22	1
EuRho/Phönix	Triglycerides, medium-chain	13020004-02	1
Gehe	Triglycerides, medium-chain	1302004	1
Henry Lamotte	Triglycerides, medium-chain	1302004-01	1
Caelo	Triglycerides, medium-chain	1302004-01	1
Alliance Healthcare	Triglycerides, medium-chain	1302004-02	1
Phönix Gotha	Triglycerides, medium-chain	1302004-02	1
Sanacorp	Triglycerides, medium-chain	130502	1
Fagron	Triglycerides, medium-chain	130502	1
Euro OTC	Triglycerides, medium-chain	1309025-01	2
Gehe	Triglycerides, medium-chain	1309025-01	1
Phoenix	Triglycerides, medium-chain	13098601	1
Caesar & Loretz GmbH	Triglycerides, medium-chain	13098601	1
Sanacorp	Triglycerides, medium-chain	13177102	1
Caelo	Triglycerides, medium-chain	13210301	1
	Triglycerides, medium-chain	13E24-N04	3
Phönix 11.11.2014	Triglycerides, medium-chain	13E2-N04	1
Fagron	Triglycerides, medium-chain	13H02-N03	1
Fagron	Triglycerides, medium-chain	13J21-N06	1
Fagron	Triglycerides, medium-chain	13j25-b11	1
Euro OTC	Triglycerides, medium-chain	1402019-01	2
Euro OTC	Triglycerides, medium-chain	141100901	1
Caesar&lorenz,Phoenix	Triglycerides, medium-chain	14219003	1
Caelo	Triglycerides, medium-chain	14K28-N01	1
Euro OTC	Triglycerides, medium-chain	150902201	1
Caelo	Triglycerides, medium-chain	15374201	1
THC Pharma	Triglycerides, medium-chain	2012050379	1
Caelo	Triglycerides, medium-chain	23082013B	1
Caelo	Triglycerides, medium-chain	3113506	3
Caelo	Triglycerides, medium-chain	60492256	1
Caelo	Triglycerides, medium-chain	AKD5	1
Fagron	Triglycerides, medium-chain	Ch120604	1
endima	Triglycerides, medium-chain	112k04	1

- 90 110 spectra from 501 *Apo-Ident* customers from a total of 26 697 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Triglycerides, medium-chain* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Triglycerides, medium-chain* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	240	0	35 304
Type B	0	550	10	93 314
Type C	0	1230	12	90 110

The substance/substance group *Triglycerides, medium-chain* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 97.5000 %)
Type B	100.0000 % (> 99.9761 %)	98.2143 % (> 97.6786 %)
Type C	100.0000 % (> 99.6586 %)	99.0338 % (> 98.7923 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31483	31483	0.00	22.56
31484	31804	1.83	19.16
31485	31804	2.11	19.26
31771	31771	0.00	21.15
31804	31804	0.00	18.59
31955	31955	0.00	19.87

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Triglycerol diisostearate**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30756-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Triglycerol diisostearate; Triglycerinum diisostearinicum; Triglyceroli diisostearas

Special notes

When selecting the *Triglycerol diisostearate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Triglycerol diisostearate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Triglycerol diis...	12018503	30756	40	20120201
Caelo	Triglycerol diis...	13081002	31219	40	1403749
Caelo	Triglycerol diis...	14118801	31338	40	20140425

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 120 spectra of 3 reference samples from the substance/substance group *Triglycerol diisostearate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 424 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 3 reference samples of the substance/substance group *Triglycerol diisostearate*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Triglycerol diisostearate	12018501	30568	40
Caelo	Triglycerol diisostearate	163311	33094	40
Caelo	Triglycerol diisostearate	180592	34073	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 132 spectra from 70 *Apo-Ident* customers from 30 batches from the substance/substance group *Triglycerol diisostearate*.
- These include spectra of independent samples from 27 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Triglycerol diisostearate	13399003	5
Caesar & Loretz GmbH	Triglycerol diisostearate	14054802	4
Caelo	Triglycerol diisostearate	14118801	8
Caelo	Triglycerol diisostearate	14193003	4
Caelo	Triglycerol diisostearate	12359401	4
Caelo	Triglycerol diisostearate	12359405	7
Caelo	Triglycerol diisostearate	13081002	5
Caelo	Triglycerol diisostearate	13236704	8
Caelo	Triglycerol diisostearate	14054802	2
gehe	Triglycerol diisostearate	14054802	1
Gehe	Triglycerol diisostearate	14118801	2
Audor Pharma	Triglycerol diisostearate	14118801	1
Caelo	Triglycerol diisostearate	14118804	1
Caelo	Triglycerol diisostearate	14193008	5
Noweda	Triglycerol diisostearate	14193008	1
Caelo	Triglycerol diisostearate	14380504	7
Noweda	Triglycerol diisostearate	14380504	1
Sanacorp	Triglycerol diisostearate	15137505	1
Caelo	Triglycerol diisostearate	15137505	10
Caelo	Triglycerol diisostearate	15283401	1
Caelo	Triglycerol diisostearate	15283402	5
Caelo	Triglycerol diisostearate	15283403	1
Caelo	Triglycerol diisostearate	15283404	7
Caesar/Noweda	Triglycerol diisostearate	15283404	1
Caelo	Triglycerol diisostearate	15283406	2
Caelo	Triglycerol diisostearate	152835405	1
Caelo	Triglycerol diisostearate	16062502	9
Phoenix	Triglycerol diisostearate	15283405	1
Caelo	Triglycerol diisostearate	16179202	1
Gehe	Triglycerol diisostearate	16062502	1
Caelo	Triglycerol diisostearate	16179205	2
Caelo	Triglycerol diisostearate	16331101	6
Caesar&Loretz GmbH	Triglycerol diisostearate	16331101	1
Caelo	Triglycerol diisostearate	12018501	2
Caelo	Triglycerol diisostearate	11162902	1
Caelo	Triglycerol diisostearate	9082013A	1
Caelo	Triglycerol diisostearate	16331107	1
Caelo	Triglycerol diisostearate	12018503	3
Caelo	Triglycerol diisostearate	12133303	5
Gehe	Triglycerol diisostearate	12359401	1
Phoenix	Triglycerol diisostearate	14054802	1
Caelo	Triglycerol diisostearate	11074801	1

- 91 220 spectra from 501 *Apo-Ident* customers from a total of 26 980 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Triglycerol diisostearate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Triglycerol diisostearate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	120	0	35 424
Type B	0	120	0	93 754
Type C	0	130	2	91 220

The substance/substance group *Triglycerol diisostearate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9705 %)	100.0000 % (> 95.0000 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6586 %)	98.4848 % (> 96.2121 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30756	31219	5.85	39.57
31219	31219	0.00	44.78
31338	31219	6.65	45.20

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Verbena oil, nature-identical
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30666-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Verbena oil, nature-identical; Oleum verbenae

Special notes

When selecting the *Verbena oil, nature-identical* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Verbena oil, nature-identical*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Verbena oil, nat...	10368808	30666	40	20101229
Caelo	Verbena oil, nat...	13141710	31194	40	1412074
Caelo	Verbena oil, nat...	14026003	31399	40	20140128
Caelo	Verbena oil, nat...	15229809	32071	40	20150714

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Verbena oil, nature-identical*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 100 spectra from 3 reference samples of the substance/substance group *Verbena oil, nature-identical*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Verbena oil, nature-identical	15229809	32071 [†]	20
Caelo	Verbena oil, nature-identical	16290903	33604	40
Caelo	Verbena oil, nature-identical	172589	33906	40

- 93 774 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1 spectra from 1 *Apo-Ident* customers from 1 batches from the substance/substance group *Verbena oil, nature-identical*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Verbena oil, nature-identical	13141704	1

- 91 351 spectra from 501 *Apo-Ident* customers from a total of 27 008 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Verbena oil, nature-identical* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Verbena oil, nature-identical* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	100	0	93 774
Type C	0	1	0	91 351

The substance/substance group *Verbena oil, nature-identical* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 94.0000 %)
Type C	100.0000 % (> 99.6675 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30666	31194	2.32	57.15
31194	31194	0.00	59.42
31399	31194	3.36	58.69
32071	31194	4.40	56.49

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Viscous base solution (NRF p. 20)**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30659-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Viscous base solution (NRF p. 20)

Special notes

When selecting the *Viscous base solution (NRF p. 20)* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Viscous base solution (NRF p. 20)*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Viscous base sol. ...	12127301	30659	40	1403145
Caelo	Viscous base sol. ...	14020401	31354	40	20140124
Caelo	Viscous base sol. ...	14109401	31658	40	AR-15-FG-007904-01
Caelo	Viscous base sol. ...	15265401	31790	40	20150825
Caelo	Viscous base sol. ...	153888	31884	40	20151124*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Viscous base sol...	153889	31886	40	20151124*
Caelo	Viscous base sol...	154121	31887	40	20151130*
Caelo	Viscous base sol...	153890	31890	40	20151124*
Caelo	Viscous base sol...	153894	31893	40	20151130*
Caelo	Viscous base sol...	153895	31896	40	20151130*
Caelo	Viscous base sol...	153896	31897	40	20151130*
Caelo	Viscous base sol...	153897	31900	40	20151130*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 480 spectra of 12 reference samples from the substance/substance group *Viscous base solution* (NRF p. 20). These samples are listed above in the section [calibration samples](#). The reference samples originate from 12 different batches.
- 35 064 spectra from a total of 828 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 2520 spectra from 63 reference samples of the substance/substance group *Viscous base solution* (NRF p. 20).
- These include spectra of independent samples from 43 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Viscous base solution (NRF p...	10160601	30161	40

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Viscous base solution (NRF p...	10232201	30167	40
Caelo	Viscous base solution (NRF p...	14020401	31354 [†]	20
Caelo	Viscous base solution (NRF p...	14109401	31658 [†]	20
Caelo	Viscous base solution (NRF p...	15265401	31790 [†]	20
Caelo	Viscous base solution (NRF p...	153888	31882	60
Caelo	Viscous base solution (NRF p...	153888	31883	60
Caelo	Viscous base solution (NRF p...	153888	31884 [†]	20
Caelo	Viscous base solution (NRF p...	153889	31885	60
Caelo	Viscous base solution (NRF p...	153889	31886 [†]	20
Caelo	Viscous base solution (NRF p...	154121	31887 [†]	20
Caelo	Viscous base solution (NRF p...	154121	31888	60
Caelo	Viscous base solution (NRF p...	153890	31889	60
Caelo	Viscous base solution (NRF p...	153890	31890 [†]	20
Caelo	Viscous base solution (NRF p...	153892	31891	60
Caelo	Viscous base solution (NRF p...	153893	31892	60
Caelo	Viscous base solution (NRF p...	153894	31893 [†]	20
Caelo	Viscous base solution (NRF p...	153894	31894	60
Caelo	Viscous base solution (NRF p...	153895	31895	60
Caelo	Viscous base solution (NRF p...	153895	31896 [†]	20
Caelo	Viscous base solution (NRF p...	153896	31897 [†]	20
Caelo	Viscous base solution (NRF p...	153896	31898	60
Caelo	Viscous base solution (NRF p...	153897	31899	60
Caelo	Viscous base solution (NRF p...	153897	31900 [†]	20
Caelo	Viscous base solution (NRF p...	161076	32287	40
Caelo	Viscous base solution (NRF p...	161089	32288	40
Caelo	Viscous base solution (NRF p...	161090	32289	40
Caelo	Viscous base solution (NRF p...	162172	32356	40
Caelo	Viscous base solution (NRF p...	162173	32357	40
Caelo	Viscous base solution (NRF p...	162174	32358	40
Caelo	Viscous base solution (NRF p...	162175	32359	40
Caelo	Viscous base solution (NRF p...	162176	32360	40
Caelo	Viscous base solution (NRF p...	162297	32554	40
Caelo	Viscous base solution (NRF p...	162298	32555	40
Caelo	Viscous base solution (NRF p...	162299	32556	40
Caelo	Viscous base solution (NRF p...	162300	32557	40
Caelo	Viscous base solution (NRF p...	162301	32558	40
Caelo	Viscous base solution (NRF p...	162650	32559	40
Caelo	Viscous base solution (NRF p...	162651	32572	40
Caelo	Viscous base solution (NRF p...	162652	32573	40
Caelo	Viscous base solution (NRF p...	162653	32574	40
Caelo	Viscous base solution (NRF p...	162655	32575	40
Caelo	Viscous base solution (NRF p...	163153	32897	40
Caelo	Viscous base solution (NRF p...	163152	32898	40
Caelo	Viscous base solution (NRF p...	170176	33146	40
Caelo	Viscous base solution (NRF p...	170177	33147	40
Caelo	Viscous base solution (NRF p...	170178	33148	40
Caelo	Viscous base solution (NRF p...	170179	33149	40
Caelo	Viscous base solution (NRF p...	170180	33150	40
Caelo	Viscous base solution (NRF p...	170239	33151	40
Caelo	Viscous base solution (NRF p...	170240	33152	40
Caelo	Viscous base solution (NRF p...	172319	33730	40
Caelo	Viscous base solution (NRF p...	171940	33731	40
Caelo	Viscous base solution (NRF p...	171941	33732	40
Caelo	Viscous base solution (NRF p...	171942	33733	40
Caelo	Viscous base solution (NRF p...	171943	33734	40
Caelo	Viscous base solution (NRF p...	171944	33735	40
Caelo	Viscous base solution (NRF p...	172892	33972	40
Caelo	Viscous base solution (NRF p...	172895	33973	40
Caelo	Viscous base solution (NRF p...	180306	34131	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Viscous base solution (NRF p...	181240	34132	40
Caelo	Viscous base solution (NRF p...	181242	34249	40
Caelo	Viscous base solution (NRF p...	181244	34250	40

- 91 354 spectra from a total of 2113 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 827 spectra from 79 *Apo-Ident* customers from 316 batches from the substance/substance group *Viscous base solution (NRF p. 20)*.
- These include spectra of independent samples from 311 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Fagron	Viscous base solution (NRF p...	1301F003	3
Apomix / 11ol	Viscous base solution (NRF p...	13229501	1
Caelo	Viscous base solution (NRF p...		1
Caelo	Viscous base solution (NRF p...	1071501	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	2012015	2
Sonnenberg-Apotheke	Viscous base solution (NRF p...	11214	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	10814	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	20516	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	40117	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	30216	1
Caelo	Viscous base solution (NRF p...	5051502	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	70514	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	70916	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	71015	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	90714	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	71216	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	91116	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	90915	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	90414	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	100314	1
Fagron	Viscous base solution (NRF p...	10031502	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	110116	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	100715	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	110516	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	110614	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	100914	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	111114	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	111115	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Sonnenberg-Apotheke	Viscous base solution (NRF p...	121016	1
Caelo	Viscous base solution (NRF p...	12127301	1
Fagron	Viscous base solution (NRF p...	1301F001	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	130416	1
Aliance am 08.04.2014	Viscous base solution (NRF p...	1306F002	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	131014	1
Fagron	Viscous base solution (NRF p...	1310167	9
Caelo	Viscous base solution (NRF p...	13229502	3
Nachprüfung	Viscous base solution (NRF p...	13158001	2
Fagron	Viscous base solution (NRF p...	13J22-T05-007729	8
Caelo	Viscous base solution (NRF p...	13157701	1
Fagron	Viscous base solution (NRF p...	1310F006	1
Caelo	Viscous base solution (NRF p...	14000801	7
Fagron	Viscous base solution (NRF p...	13J22-T05-007720	1
Caelo	Viscous base solution (NRF p...	14000901	2
Caelo	Viscous base solution (NRF p...	14020301	5
Caelo	Viscous base solution (NRF p...	14020401	10
VDL 24.05.14	Viscous base solution (NRF p...	14020401	2
Caelo	Viscous base solution (NRF p...	14020402	6
Fagron	Viscous base solution (NRF p...	14020502	5
Gutenberg-Apotheke	Viscous base solution (NRF p...	14020502	1
GEHE am 27.05.2014	Viscous base solution (NRF p...	14020502	1
Caelo	Viscous base solution (NRF p...	14020701	1
Fagron	Viscous base solution (NRF p...	1402F001	25
Caelo	Viscous base solution (NRF p...	14058901	3
Caelo	Viscous base solution (NRF p...	14059001	3
Fagron	Viscous base solution (NRF p...	14059001	2
Caelo	Viscous base solution (NRF p...	14059101	7
Gutenberg-Apotheke	Viscous base solution (NRF p...	14059101	3
Caelo	Viscous base solution (NRF p...	14059301	1
Caelo	Viscous base solution (NRF p...	14109401	6
VDL 17.10.14	Viscous base solution (NRF p...	14109401	2
Caesar und Loretz	Viscous base solution (NRF p...	14109401	1
Caelo	Viscous base solution (NRF p...	14109402	3
Sanacorp	Viscous base solution (NRF p...	14109402	1
Caelo	Viscous base solution (NRF p...	14109501	5
Caelo	Viscous base solution (NRF p...	14109801	3
Caelo	Viscous base solution (NRF p...	14109701	1
Caelo	Viscous base solution (NRF p...	14109901	11
Alliance Health.18.09.2...	Viscous base solution (NRF p...	14109901	1
Phoenix 08.07.2014	Viscous base solution (NRF p...	14110001	1
Caelo	Viscous base solution (NRF p...	14110001	2
Caelo	Viscous base solution (NRF p...	14187201	9
Caelo	Viscous base solution (NRF p...	1418730	2
Caelo	Viscous base solution (NRF p...	14187301	8
Caelo	Viscous base solution (NRF p...	14187501	1
Caelo	Viscous base solution (NRF p...	14187601	5
Caelo	Viscous base solution (NRF p...	14187801	12
Gehe 19.12.14	Viscous base solution (NRF p...	14187601	1
Caelo	Viscous base solution (NRF p...	14187701	1
Caelo	Viscous base solution (NRF p...	14187901	3
Caelo	Viscous base solution (NRF p...	14188001	5
Caelo	Viscous base solution (NRF p...	14266101	1
Caelo	Viscous base solution (NRF p...	14266501	7
Caelo	Viscous base solution (NRF p...	14266601	2
Caelo	Viscous base solution (NRF p...	14354304	3
Caelo	Viscous base solution (NRF p...	14354401	5
Caelo	Viscous base solution (NRF p...	14368401	5
Caelo	Viscous base solution (NRF p...	14354501	1
Caelo	Viscous base solution (NRF p...	14368601	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Viscous base solution (NRF p...	14368701	3
Fagron	Viscous base solution (NRF p...	14K24-T03-031283	2
Caelo	Viscous base solution (NRF p...	14368801	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	150616	1
Fagron	Viscous base solution (NRF p...	1509255	2
Caelo	Viscous base solution (NRF p...	15093301	10
Caelo	Viscous base solution (NRF p...	15093401	4
Caelo	Viscous base solution (NRF p...	15093501	2
Caelo	Viscous base solution (NRF p...	15093601	4
Caelo	Viscous base solution (NRF p...	15093602	4
Caelo	Viscous base solution (NRF p...	15093701	3
Caelo	Viscous base solution (NRF p...	15093801	2
Caelo	Viscous base solution (NRF p...	15105701	4
Caelo	Viscous base solution (NRF p...	15141601	2
Euro OTC	Viscous base solution (NRF p...	15141601	1
Caelo	Viscous base solution (NRF p...	15141901	2
Caelo	Viscous base solution (NRF p...	15180101	3
Caelo	Viscous base solution (NRF p...	15142001	1
Caelo	Viscous base solution (NRF p...	15196601	1
Caelo	Viscous base solution (NRF p...	15196701	2
Caelo	Viscous base solution (NRF p...	15264801	7
Caelo	Viscous base solution (NRF p...	152364901	1
Caelo	Viscous base solution (NRF p...	15264901	1
Caelo	Viscous base solution (NRF p...	15265001	1
Caelo	Viscous base solution (NRF p...	15265201	4
Caelo	Viscous base solution (NRF p...	15265101	2
Caelo	Viscous base solution (NRF p...	15265301	9
Sanacorp	Viscous base solution (NRF p...	15265201	1
Caelo	Viscous base solution (NRF p...	15265401	8
Fishar	Viscous base solution (NRF p...	15265401	1
Caelo	Viscous base solution (NRF p...	15265402	1
Caelo	Viscous base solution (NRF p...	15265501	6
Caelo	Viscous base solution (NRF p...	15265601	2
Sanacorp	Viscous base solution (NRF p...	15265601	1
Caelo	Viscous base solution (NRF p...	15265701	9
Caesar&Loretz GmbH	Viscous base solution (NRF p...	15265701	1
Caelo	Viscous base solution (NRF p...	15265801	3
Caelo	Viscous base solution (NRF p...	15265901	5
Caelo	Viscous base solution (NRF p...	15389401	3
Caelo	Viscous base solution (NRF p...	15389501	2
Caelo	Viscous base solution (NRF p...	15389701	7
Sonnenberg-Apotheke	Viscous base solution (NRF p...	160415	1
Fagron	Viscous base solution (NRF p...	1604682	3
Sonnenberg-Apotheke	Viscous base solution (NRF p...	160316	1
Caelo	Viscous base solution (NRF p...	15412101	2
Defaktur Bahnhof Apotheke	Viscous base solution (NRF p...	160630008	1
Defaktur Bahnhof Apotheke	Viscous base solution (NRF p...	160630007	1
Defaktur Bahnhof Apotheke	Viscous base solution (NRF p...	160718007	1
Defaktur Bahnhof Apotheke	Viscous base solution (NRF p...	160719002	1
Defaktur Bahnhof Apotheke	Viscous base solution (NRF p...	160805001	1
Defaktur Bahnhof Apotheke	Viscous base solution (NRF p...	160805002	1
Bahnhof Apotheke	Viscous base solution (NRF p...	160818002	1
Bahnhof Apotheke	Viscous base solution (NRF p...	160818001	1
Bahnhof-Apotheke	Viscous base solution (NRF p...	160825001	1
Bahnhof-Apotheke	Viscous base solution (NRF p...	160823001	1
Caesar&Loretz GmbH	Viscous base solution (NRF p...	16107601	3
Bahnhof-Apotheke	Viscous base solution (NRF p...	160825002	1
Caelo	Viscous base solution (NRF p...	16107601	4
Caelo	Viscous base solution (NRF p...	16108001	3
Caelo	Viscous base solution (NRF p...	16108201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Viscous base solution (NRF p...	16108101	1
Caelo	Viscous base solution (NRF p...	16108301	3
Caelo	Viscous base solution (NRF p...	16108401	2
Caelo	Viscous base solution (NRF p...	16108501	3
Caelo	Viscous base solution (NRF p...	16108601	4
sanacorp, 32,14EUR, 15.09...	Viscous base solution (NRF p...	16108601	1
Caelo	Viscous base solution (NRF p...	16108901	8
Jenne	Viscous base solution (NRF p...	16108901	1
Caesar&Loretz GmbH	Viscous base solution (NRF p...	16108901	1
Caelo	Viscous base solution (NRF p...	16108701	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	161215	1
Caelo	Viscous base solution (NRF p...	16109001	1
Caelo	Viscous base solution (NRF p...	16159401	1
Caelo	Viscous base solution (NRF p...	16217301	5
Caesar&Loretz GmbH	Viscous base solution (NRF p...	16217301	1
Caelo	Viscous base solution (NRF p...	16217601	4
Apomix/Phönix	Viscous base solution (NRF p...	16217601	1
Caelo	Viscous base solution (NRF p...	16217602	10
Caelo	Viscous base solution (NRF p...	16230001	2
Caesar&Loretz GmbH	Viscous base solution (NRF p...	16229701	1
Caelo	Viscous base solution (NRF p...	16230101	1
Caelo	Viscous base solution (NRF p...	16265001	2
Caelo	Viscous base solution (NRF p...	16315201	4
Caelo	Viscous base solution (NRF p...	16265101	1
Sanacorp, 22.02.17, 32,...	Viscous base solution (NRF p...	16265101	1
Caelo	Viscous base solution (NRF p...	16315301	5
Caelo	Viscous base solution (NRF p...	17017801	2
Caelo	Viscous base solution (NRF p...	17023901	1
Fagron	Viscous base solution (NRF p...	1705701	2
gehe	Viscous base solution (NRF p...	17097801	1
Caelo	Viscous base solution (NRF p...	17097901	1
Caelo	Viscous base solution (NRF p...	17098101	1
Caelo	Viscous base solution (NRF p...	17098102	4
Caelo	Viscous base solution (NRF p...	17194201	1
gehe	Viscous base solution (NRF p...	17194401	2
Caelo	Viscous base solution (NRF p...	17193502	1
gehe	Viscous base solution (NRF p...	17194101	1
Caelo	Viscous base solution (NRF p...	18071502	1
Caelo	Viscous base solution (NRF p...	17231201	1
Caelo	Viscous base solution (NRF p...	17194401	1
Caelo	Viscous base solution (NRF p...	17138401	1
Caelo	Viscous base solution (NRF p...	17231301	1
Caelo	Viscous base solution (NRF p...	17138801	1
Euro OTC	Viscous base solution (NRF p...	17232001	1
Caelo	Viscous base solution (NRF p...	20081501	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	211216	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	220616	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	220216	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	230315	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	250215	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	260115	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	250716	1
Fagron	Viscous base solution (NRF p...	28111402	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	280916	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	290716	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	231116	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	270515	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	290716-2	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	261016	1
Sonnenberg Apotheke	Viscous base solution (NRF p...	300615	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sonnenberg-Apotheke	Viscous base solution (NRF p...	60114	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	91213	1
Caelo	Viscous base solution (NRF p...	3x14058901	1
Sonnenberg-Apotheke	Viscous base solution (NRF p...	100214	1
Fagron	Viscous base solution (NRF p...	1111F001	1
Caelo	Viscous base solution (NRF p...	11140501	1
Klenk	Viscous base solution (NRF p...	1111F003	1
Fagron	Viscous base solution (NRF p...	1111F003	1
Anzag	Viscous base solution (NRF p...	12018901	1
Noweda	Viscous base solution (NRF p...	11176201	1
Caelo	Viscous base solution (NRF p...	11176401	1
Caelo	Viscous base solution (NRF p...	12019001	1
Caelo	Viscous base solution (NRF p...	12019002	2
Noweda	Viscous base solution (NRF p...	12019401	1
Gehe	Viscous base solution (NRF p...	12019201	1
Fagron	Viscous base solution (NRF p...	1203C004	2
Fagron	Viscous base solution (NRF p...	1206F003	3
Noweda	Viscous base solution (NRF p...	12093601	2
Caelo	Viscous base solution (NRF p...	12093701	1
Caelo	Viscous base solution (NRF p...	12094101	2
Caelo	Viscous base solution (NRF p...	12093901	1
Caelo	Viscous base solution (NRF p...	12093801	1
Anzag 05.09.2012	Viscous base solution (NRF p...	12127302	4
Sanacorp	Viscous base solution (NRF p...	12127401	2
Caelo	Viscous base solution (NRF p...	12127501	1
Caelo	Viscous base solution (NRF p...	12127601	2
Anzag am 28.01.2013	Viscous base solution (NRF p...	12127601	2
Fagron	Viscous base solution (NRF p...	1214509	1
Caelo	Viscous base solution (NRF p...	12231101	5
Caelo	Viscous base solution (NRF p...	12231002	1
Sanacorp	Viscous base solution (NRF p...	12231301	4
Caelo	Viscous base solution (NRF p...	12231201	1
Caelo	Viscous base solution (NRF p...	12249801	1
Caelo	Viscous base solution (NRF p...	12249901	3
Caelo	Viscous base solution (NRF p...	1224991	1
Phönix, 27.12.2012	Viscous base solution (NRF p...	12250101	3
Caelo	Viscous base solution (NRF p...	12309001	2
Caelo	Viscous base solution (NRF p...	12309101	11
Caelo	Viscous base solution (NRF p...	12309201	6
Caelo	Viscous base solution (NRF p...	12309202	3
Phoenix 11.03.2013	Viscous base solution (NRF p...	12309202	5
Apomix / 11	Viscous base solution (NRF p...	12309201	1
Phönix 04.04.13	Viscous base solution (NRF p...	13006001	6
Caelo	Viscous base solution (NRF p...	13006001	3
Caelo	Viscous base solution (NRF p...	13006201	4
Anzag/Hedinger	Viscous base solution (NRF p...	13006301	5
Caelo	Viscous base solution (NRF p...	13006401	6
Caelo	Viscous base solution (NRF p...	13006501	2
Fagron	Viscous base solution (NRF p...	1301673	4
Phönix 11.11.2014	Viscous base solution (NRF p...	1301F001	1
Fagron	Viscous base solution (NRF p...	1301F002	3
Fagron	Viscous base solution (NRF p...	1302F002	18
Caelo	Viscous base solution (NRF p...	13096001	2
Post Aptotheke	Viscous base solution (NRF p...	1306F002	2
Caelo	Viscous base solution (NRF p...	13071302	1
Caelo	Viscous base solution (NRF p...	13096101	7
Herstellung	Viscous base solution (NRF p...	131018001	1
Herstellung	Viscous base solution (NRF p...	131018004	1
Herstellung	Viscous base solution (NRF p...	131025001	1
Herstellung	Viscous base solution (NRF p...	131030004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Herstellung	Viscous base solution (NRF p...	131115001	1
Herstellung	Viscous base solution (NRF p...	131030005	1
Herstellung	Viscous base solution (NRF p...	131126002	2
Herstellung	Viscous base solution (NRF p...	131206001	1
Herstellung	Viscous base solution (NRF p...	131203012	1
Herstellung	Viscous base solution (NRF p...	131211001	1
Herstellung	Viscous base solution (NRF p...	131211002	1
Herstellung	Viscous base solution (NRF p...	131218006	1
Herstellung	Viscous base solution (NRF p...	131220001	1
Caelo	Viscous base solution (NRF p...	13150801	3
Post-Apotheke	Viscous base solution (NRF p...	13150801	3
Caelo	Viscous base solution (NRF p...	13151101	6
Caelo	Viscous base solution (NRF p...	13150802	2
Caelo	Viscous base solution (NRF p...	13151201	4
Caelo	Viscous base solution (NRF p...	13158001	8
Caelo	Viscous base solution (NRF p...	13158002	3
Caelo	Viscous base solution (NRF p...	1315801	1
Caelo	Viscous base solution (NRF p...	13158101	2
Caesar & Loretz GmbH	Viscous base solution (NRF p...	13158201	1
Phönix 11.11.2013	Viscous base solution (NRF p...	13158101	1
Anzag	Viscous base solution (NRF p...	13158201	2
Caelo	Viscous base solution (NRF p...	13158201	2
Caelo	Viscous base solution (NRF p...	13220101	1
Caelo	Viscous base solution (NRF p...	13228901	8
Caelo	Viscous base solution (NRF p...	13229001	4
Caelo	Viscous base solution (NRF p...	13229101	1
Caelo	Viscous base solution (NRF p...	13229201	2
Caelo	Viscous base solution (NRF p...	13229501	4
Alliance Health 21.02.2...	Viscous base solution (NRF p...	13229502	2
Caelo	Viscous base solution (NRF p...	13246101	1
Caelo	Viscous base solution (NRF p...	13246401	5
Caelo	Viscous base solution (NRF p...	13246501	6
Caelo	Viscous base solution (NRF p...	13246701	5
Caelo	Viscous base solution (NRF p...	13246901	1
Herstellung	Viscous base solution (NRF p...	140108002	1
Herstellung	Viscous base solution (NRF p...	140117001	1
Herstellung	Viscous base solution (NRF p...	140108003	1
Caelo	Viscous base solution (NRF p...	14020601	3
Bombastus	Viscous base solution (NRF p...	15265501	1
Caelo	Viscous base solution (NRF p...	16107901	1
Herstellung	Viscous base solution (NRF p...	20030513-1	1
Herstellung	Viscous base solution (NRF p...	20010213-1	1
Herstellung	Viscous base solution (NRF p...	20040113-2	1
Herstellung	Viscous base solution (NRF p...	20040113-1	1
Herstellung	Viscous base solution (NRF p...	20041013-1	1
Herstellung	Viscous base solution (NRF p...	20050613-1	1
Herstellung	Viscous base solution (NRF p...	20050713-1	1
Herstellung	Viscous base solution (NRF p...	20060913-1	1
Herstellung	Viscous base solution (NRF p...	20050713-2	1
Herstellung	Viscous base solution (NRF p...	20060913-2	1
Herstellung	Viscous base solution (NRF p...	20080213-1	1
Herstellung	Viscous base solution (NRF p...	20090113-1	2
Herstellung	Viscous base solution (NRF p...	20080313-1	1
Herstellung	Viscous base solution (NRF p...	20100513-1	1
Herstellung	Viscous base solution (NRF p...	20100613-1	1
Herstellung	Viscous base solution (NRF p...	20120713-1	1
Herstellung	Viscous base solution (NRF p...	20150213-1	1
Herstellung	Viscous base solution (NRF p...	20130613-1	1
Herstellung	Viscous base solution (NRF p...	20130313-1	1
Herstellung	Viscous base solution (NRF p...	20150813-4	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Herstellung	Viscous base solution (NRF p...	20160513-2	1
Herstellung	Viscous base solution (NRF p...	20140613-1	1
Herstellung	Viscous base solution (NRF p...	20150813-3	1
Herstellung	Viscous base solution (NRF p...	20190413-1	1
Herstellung	Viscous base solution (NRF p...	20200913-1	1
Herstellung	Viscous base solution (NRF p...	20190413-2	1
Herstellung	Viscous base solution (NRF p...	20210313-2	1
Herstellung	Viscous base solution (NRF p...	20210313-1	1
Herstellung	Viscous base solution (NRF p...	20220213-1	1
Herstellung	Viscous base solution (NRF p...	20220513-2	1
Herstellung	Viscous base solution (NRF p...	20220213-2	1
Caelo	Viscous base solution (NRF p...	20250113-1	1
Herstellung	Viscous base solution (NRF p...	20260313-1	1
Herstellung	Viscous base solution (NRF p...	20260313-2	1
Herstellung	Viscous base solution (NRF p...	20260713-2	1
Herstellung	Viscous base solution (NRF p...	20260713-1	1
Herstellung	Viscous base solution (NRF p...	20270913-2	1
Herstellung	Viscous base solution (NRF p...	20280513-1	1
Herstellung	Viscous base solution (NRF p...	20280613-1	1
Herstellung	Viscous base solution (NRF p...	20300713-2	1
Herstellung	Viscous base solution (NRF p...	20280613-2	1
Herstellung	Viscous base solution (NRF p...	2030913-1	1
Herstellung	Viscous base solution (NRF p...	2150313-1	1
Fagron	Viscous base solution (NRF p...	22041304	2
Herstellung	Viscous base solution (NRF p...	240114	1
Caelo	Viscous base solution (NRF p...	3013A-03142	1
Fagron	Viscous base solution (NRF p...	31101205	1
Anzag, 18.07.2012	Viscous base solution (NRF p...		1

- 90 525 spectra from 501 *Apo-Ident* customers from a total of 26 698 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Viscous base solution (NRF p. 20)* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Viscous base solution (NRF p. 20)* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	480	0	35 064
Type B	0	2510	10	91 354
Type C	0	793	34	90 525

The substance/substance group *Viscous base solution (NRF p. 20)* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.7500 %)
Type B	100.0000 % (> 99.9761 %)	99.6032 % (> 99.4841 %)
Type C	100.0000 % (> 99.6586 %)	95.8888 % (> 95.5260 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30659	30659	0.00	26.18
31354	31900	1.60	23.58
31658	31658	0.00	27.38
31790	31893	4.16	27.52
31884	31884	0.00	27.56
31886	31886	0.00	27.17
31887	31887	0.00	28.15
31890	31890	0.00	27.63
31893	31893	0.00	28.70
31896	31896	0.00	27.08
31897	31897	0.00	26.72
31900	31900	0.00	27.73

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Water and aqueous solutions
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30869-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Water and aqueous solutions; Aluminium acetotartrate solution; Aqua; Chlorhexidine digluconate solution 20%; Chlorhexidini digluconatis solutio 20%; Liquor aluminii acetico-tartarici; Solutio chlorhexidini gluconici 20%; Water

Special notes

When selecting the *Water and aqueous solutions* substance/substance group, the following information is displayed to the user:

- Advice for the identification of water: Preserved, purified and distilled water can not be conclusively distinguished by spectral analysis. An additional analysis of identity may be required.

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Water and aqueous solutions*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Aluminium acetot...	12233502	30869	40	20120810
Caelo	Water	13176101	31057	40	20130627

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Chlorhexidine di...	13H19-B01	31240	30	20130906
Caelo	Aluminium acetot...	13378704	31342	40	AR-14-FG-012698-01
Caelo	Chlorhexidine di...	14011009	31463	40	AR-15-FG-010164-03
Caelo	Water	14253301	31624	40	20140923
Fagron	Water	14K05-T02-016775	31676	40	20150225
Caelo	Chlorhexidine di...	15131504	31705	40	20150424
Caelo	Water	15142501	31718	40	1510212
Euro OTC	Chlorhexidine di...	1508022	31735	40	20150916*
Fagron	Water	15D01-C20-16500	31809	40	20151105
Fagron	Chlorhexidine di...	15I03-B02-315213	31829	30	20151027
Euro OTC	Chlorhexidine di...	1511039	31846	40	20160114*
Caelo	Aluminium acetot...	15299901	31959	40	20151001
Caelo	Chlorhexidine di...	15292508	32308	40	20151109
Caelo	Water	161975	32354	30	20160725*
Euro OTC	Chlorhexidine di...	1511033-01	32732	40	20160114
Euro OTC	Chlorhexidine di...	1610021	32831	40	20161121*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 690 spectra of 18 reference samples from the substance/substance group *Water and aqueous solutions*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 18 different batches.
- 34 854 spectra from a total of 822 batches from further 209 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 1575 spectra from 45 reference samples of the substance/substance group *Water and aqueous solutions*.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

- These include spectra of independent samples from 35 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Chlorhexidine digluconate so...	82369169	30005	20
Fagron	Chlorhexidine digluconate so...	11C16-N07	30246	40
Fagron	Chlorhexidine digluconate so...	10H24-N08	30268	40
Caelo	Chlorhexidine digluconate so...	14011009	31463 [†]	20
Caelo	Water	14253301	31624 [†]	20
Fagron	Water	14K05-T02-016775	31676 [†]	20
Caelo	Chlorhexidine digluconate so...	15131504	31705 [†]	20
Caelo	Water	15142501	31718 [†]	20
Euro OTC	Chlorhexidine digluconate so...	1508022	31735 [†]	20
Fagron	Water	15D01-C20-16500	31809 [†]	20
Fagron	Chlorhexidine digluconate so...	15I03-B02-315213	31829 [†]	15
Euro OTC	Chlorhexidine digluconate so...	1511039	31846 [†]	20
Caelo	Aluminium acetotartrate solu...	15299901	31959 [†]	20
Caelo	Water	153714	32080	40
Caelo	Water	153945	32082	40
Caelo	Water	162119	32353	40
Caelo	Water	162230	32355	40
Euro OTC	Chlorhexidine digluconate so...	1610022	32832	40
Caelo	Water	163133	32893	40
Caelo	Water	162554	32894	40
Hedinger	Aluminium acetotartrate solu...	020160B	32921	40
Fagron	Water	16K25-C20-17409	32963	40
Caelo	Chlorhexidine digluconate so...	163088	33089	40
Caelo	Water	170020	33178	40
Caelo	Water	170324	33179	40
Euro OTC	Chlorhexidine digluconate so...	1701037	33230	40
Caelo	Aluminium acetotartrate solu...	16004505	33490	40
Caelo	Water	171731	33575	40
Caelo	Water	171451	33576	40
Caelo	Water	171773	33577	40
Caelo	Chlorhexidine digluconate so...	172101	33693	40
Caelo	Water	172534	33738	40
Caelo	Water	172535	33739	40
Caelo	Water	172532	33740	40
Euro OTC	Chlorhexidine digluconate so...	1708006	33843	40
Euro OTC	Chlorhexidine digluconate so...	1712002	33927	40
Caelo	Water	170988	33976	40
Caelo	Water	172958	33979	40
Caelo	Water	180004	34135	40
Caelo	Water	181355	34138	40
Euro OTC	Chlorhexidine digluconate so...	1807014	34191	40
Caelo	Water	181951	34240	40
Caelo	Water	181352	34241	40
Caelo	Water	182287	34242	40
Caelo	Water	182399	34245	40

- 92 299 spectra from a total of 2122 batches from further 434 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 2483 spectra from 386 *Apo-Ident* customers from 677 batches from the substance/substance group *Water and aqueous solutions*.
- These include spectra of independent samples from 660 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Hedinger/Phoenix	Aluminium acetotartrate solu...	20020	1
Franz GmbH	Aluminium acetotartrate solu...	20152B	2
Hedinger/Noweda	Aluminium acetotartrate solu...	200161B	1
Hedinger/Gehe	Aluminium acetotartrate solu...	20153B	4
Hedinger	Aluminium acetotartrate solu...	20154B	11
Phönix/Hedinger	Aluminium acetotartrate solu...	20155B	3
Hedinger/Noweda	Aluminium acetotartrate solu...	20155B	1
Hedinger	Aluminium acetotartrate solu...	20155B	10
Gehe	Aluminium acetotartrate solu...	20155B	1
Caelo	Aluminium acetotartrate solu...	20155B	5
Sanacorp	Aluminium acetotartrate solu...	20155B	1
Franz GmbH	Aluminium acetotartrate solu...	20155B	3
Hedinger / Alliance Hea...	Aluminium acetotartrate solu...	20155B	1
Hedinger / Sanacorp	Aluminium acetotartrate solu...	20155B	1
Hedinger über Phoenix	Aluminium acetotartrate solu...	20155B	1
Caelo	Aluminium acetotartrate solu...	201563	1
Hedinger/Noweda	Aluminium acetotartrate solu...	201568	1
Hedinger	Aluminium acetotartrate solu...	201568	1
Hedinger	Aluminium acetotartrate solu...	20156B	13
Gehe Kassel	Aluminium acetotartrate solu...	201568B	1
Raben-Apotheke Mildenau	Aluminium acetotartrate solu...	20156B	5
Hedinger/ Noweda	Aluminium acetotartrate solu...	20156B	1
Caelo	Aluminium acetotartrate solu...	20156B	4
Hedinger / Gehe	Aluminium acetotartrate solu...	20156B	4
Hedinger/Jenne 16.08.14	Aluminium acetotartrate solu...	20156B	1
hedinger/Gehe	Aluminium acetotartrate solu...	20156B	1
Hedinger/Noweda	Aluminium acetotartrate solu...	20156B	3
Hediger/Phönix	Aluminium acetotartrate solu...	20156B	1
Hedinger / Sanacorp	Aluminium acetotartrate solu...	20156B	1
Hedinger/Phönix	Aluminium acetotartrate solu...	20156B	1
hedinger	Aluminium acetotartrate solu...	20156B	3
Noweda/Hedinger	Aluminium acetotartrate solu...	20156B	1
sanacorp/Hedinger	Aluminium acetotartrate solu...	20157B	1
Hedinger/Phoenix	Aluminium acetotartrate solu...	201578	1
Hedinger	Aluminium acetotartrate solu...	20157B	8
Hedinger/Gehe	Aluminium acetotartrate solu...	20157B	1
Hedinger / Gehe	Aluminium acetotartrate solu...	20157B	1
hedinger/phönix	Aluminium acetotartrate solu...	20157B	1
Hedinger/ Noweda	Aluminium acetotartrate solu...	20157B	1
Hedinger über Phoenix	Aluminium acetotartrate solu...	20157B	1
hedinger	Aluminium acetotartrate solu...	20157B	2
GEHE 08.05	Aluminium acetotartrate solu...	20157B	1
Sanacorp/Hedinger	Aluminium acetotartrate solu...	20157B	1
Hedinger/Phönix	Aluminium acetotartrate solu...	20157B	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
hedinger	Aluminium acetotartrate solu...	20158B	3
Spangro/Hedinger	Aluminium acetotartrate solu...	20158B	1
Hedinger/Anzag	Aluminium acetotartrate solu...	20157B	1
Hedinger/Phoenix	Aluminium acetotartrate solu...	20158B	2
Hedinger/Sanacorp	Aluminium acetotartrate solu...	20158B	1
Hedinger/ Sanacorp WE...	Aluminium acetotartrate solu...	20158B	1
Hedinger/Gehe	Aluminium acetotartrate solu...	20158B	2
Hedinger/Noweda	Aluminium acetotartrate solu...	20158B	2
Hedinger / Noweda	Aluminium acetotartrate solu...	201598	1
Hedinger	Aluminium acetotartrate solu...	20158B	3
hedinger	Aluminium acetotartrate solu...	201598	1
Hedinger/Phoenix	Aluminium acetotartrate solu...	20159B	3
Hedinger/ Phönix	Aluminium acetotartrate solu...	20159B	1
Caelo	Aluminium acetotartrate solu...	20159B	2
Hedinger/Alliance Healt...	Aluminium acetotartrate solu...	20159B	1
Hedinger/Phönix	Aluminium acetotartrate solu...	20159B	4
Hedinger/ Noweda	Aluminium acetotartrate solu...	20159B	1
HEDINGER/NOWEDA	Aluminium acetotartrate solu...	20159B	1
Hedinger/Sana	Aluminium acetotartrate solu...	20159B	1
heidinger/phönix	Aluminium acetotartrate solu...	20159B	2
AHD/Hedinger	Aluminium acetotartrate solu...	20159B	1
Hedinger	Aluminium acetotartrate solu...	20159B	4
Hedinger/Gehe	Aluminium acetotartrate solu...	20159B	1
Hedinger/Anzag	Aluminium acetotartrate solu...	20159B	1
Hedinger/ Sanacorp	Aluminium acetotartrate solu...	20159B	1
Hedinger	Aluminium acetotartrate solu...	201608	1
Hedinger	Aluminium acetotartrate solu...	20160B	10
Hedinger/Phoenix	Aluminium acetotartrate solu...	20160B	3
Caelo	Aluminium acetotartrate solu...	20160B	3
WE:21.03.2016 Sanacor...	Aluminium acetotartrate solu...	20160B	1
Kehr/gepepharm	Aluminium acetotartrate solu...	20160B	1
Hedinger / Anzag	Aluminium acetotartrate solu...	20160B	1
Gehe	Aluminium acetotartrate solu...	20160B	2
Hedinger/Gehe	Aluminium acetotartrate solu...	20160B	1
Haedinger	Aluminium acetotartrate solu...	20160B	1
Hedinger/Noweda	Aluminium acetotartrate solu...	20160B	2
hedinger	Aluminium acetotartrate solu...	20160B	1
Hedinger/Phönix	Aluminium acetotartrate solu...	20160B	1
Hedinger GmbH	Aluminium acetotartrate solu...	20160B	1
L 20/Hedinger	Aluminium acetotartrate solu...	201618	1
Hedinger / Noweda	Aluminium acetotartrate solu...	20160B	1
Caelo	Aluminium acetotartrate solu...	20161B	6
Hedinger / Noweda	Aluminium acetotartrate solu...	20161B	1
Hedinger	Aluminium acetotartrate solu...	20161B	7
Franz GmbH / Phoenix	Aluminium acetotartrate solu...	20161B	1
Jenne	Aluminium acetotartrate solu...	20161B	1
Hedinger/ Sanacorp	Aluminium acetotartrate solu...	20161B	2
Hedinger/Fiebig	Aluminium acetotartrate solu...	20161B	1
Hedinger/ Sanacorp	Aluminium acetotartrate solu...	20161B	1
Aug. Hedinger GmbH & Co...	Aluminium acetotartrate solu...	20161B	1
Hedinger/Gehe	Aluminium acetotartrate solu...	20161B	3
Hedinger GmbH	Aluminium acetotartrate solu...	20161B	1
Phoenix WE 11.03.17 EK...	Aluminium acetotartrate solu...	20161B	1
GEHE	Aluminium acetotartrate solu...	20161B	1
Hedinger/Anzag	Aluminium acetotartrate solu...	20161B	1
Hedinger/VDL	Aluminium acetotartrate solu...	20161B	1
A. Hedinger/Noweda	Aluminium acetotartrate solu...	201628	1
Hedinger/Fiebig	Aluminium acetotartrate solu...	20162B	2
Franz GmbH, Phoenix	Aluminium acetotartrate solu...	20162B	2
AHD	Aluminium acetotartrate solu...	20162B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Aluminium acetotartrate solu. . .	20162B	1
Hedinger/ Gehe	Aluminium acetotartrate solu. . .	20162B	1
Hedinger	Aluminium acetotartrate solu. . .	20163b	1
Hedinger/Sanacorp	Aluminium acetotartrate solu. . .	20163B	1
Hedinger	Aluminium acetotartrate solu. . .	2016b	1
Caelo	Aluminium acetotartrate solu. . .	2016B	1
Hedinger	Aluminium acetotartrate solu. . .	5081403	1
Caelo	Aluminium acetotartrate solu. . .	6031415	1
athenstaedt GmbH / Noweda	Aluminium acetotartrate solu. . .	110024	1
athenstaedt	Aluminium acetotartrate solu. . .	110024	1
Caelo	Aluminium acetotartrate solu. . .	110029	1
Caelo	Aluminium acetotartrate solu. . .	13261302	2
Caelo	Aluminium acetotartrate solu. . .	12318703	1
Caelo	Aluminium acetotartrate solu. . .	13261303	1
Noweda	Aluminium acetotartrate solu. . .	13261304	1
Caelo	Aluminium acetotartrate solu. . .	13261304	3
Phönix Gotha	Aluminium acetotartrate solu. . .	13261304	1
Caelo	Aluminium acetotartrate solu. . .	13378701	2
Euro OTC	Aluminium acetotartrate solu. . .	13378701	1
Caelo	Aluminium acetotartrate solu. . .	13378703	1
Caelo	Aluminium acetotartrate solu. . .	13378704	1
Caelo	Aluminium acetotartrate solu. . .	13378706	3
Caelo	Aluminium acetotartrate solu. . .	13378708	2
Caelo	Aluminium acetotartrate solu. . .	14071503	1
Caelo	Aluminium acetotartrate solu. . .	13378709	3
Caelo	Aluminium acetotartrate solu. . .	1425010	1
Caelo	Aluminium acetotartrate solu. . .	1425006	1
Caelo	Aluminium acetotartrate solu. . .	14252004	1
Caelo	Aluminium acetotartrate solu. . .	14252007	1
Caesar u Loretz	Aluminium acetotartrate solu. . .	14252005	1
Caelo	Aluminium acetotartrate solu. . .	14252008	2
Caelo	Aluminium acetotartrate solu. . .	14252009	1
Caelo	Aluminium acetotartrate solu. . .	14252010	6
Fagron	Aluminium acetotartrate solu. . .	14252010	1
Caesar&Loretz GmbH/Anzag	Aluminium acetotartrate solu. . .	14252010	1
Caelo	Aluminium acetotartrate solu. . .	15045403	4
Caelo	Aluminium acetotartrate solu. . .	1504506	1
Caelo	Aluminium acetotartrate solu. . .	15045405	5
Caelo	Aluminium acetotartrate solu. . .	15045406	10
Caelo	Aluminium acetotartrate solu. . .	15299901	11
Caelo/Anzag	Aluminium acetotartrate solu. . .	15299901	1
Caelo	Aluminium acetotartrate solu. . .	15299902	4
Caelo	Aluminium acetotartrate solu. . .	15299904	3
Caelo	Aluminium acetotartrate solu. . .	15299909	4
Caelo/AHD	Aluminium acetotartrate solu. . .	15299909	1
Bombastus	Aluminium acetotartrate solu. . .	160029	1
athenstaedt/AllianceHC	Aluminium acetotartrate solu. . .	160029	1
Athenstaedt	Aluminium acetotartrate solu. . .	160029	1
Athenstaedt	Aluminium acetotartrate solu. . .	160036	2
athenstaedt / Gehe	Aluminium acetotartrate solu. . .	160036	1
Caelo	Aluminium acetotartrate solu. . .	16004501	2
Caelo	Aluminium acetotartrate solu. . .	16004504	2
Caesar & Loretz GmbH	Aluminium acetotartrate solu. . .	16004504	1
Caelo	Aluminium acetotartrate solu. . .	16004505	8
Caesar & Loretz GmbH	Aluminium acetotartrate solu. . .	16004505	1
Caelo	Aluminium acetotartrate solu. . .	16004512	4
Caelo	Aluminium acetotartrate solu. . .	16342001	2
Caelo	Aluminium acetotartrate solu. . .	16342006	2
Hedinger/Phönix	Aluminium acetotartrate solu. . .	(020)020	1
Hedinger/ Noweda	Aluminium acetotartrate solu. . .	20149B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Aluminium acetotartrate solu...	20150B	2
Hedinger, Heiligenwiese...	Aluminium acetotartrate solu...	20150B	1
Kehr	Aluminium acetotartrate solu...	20150B	1
Hedinger	Aluminium acetotartrate solu...	20149B	1
Caelo	Aluminium acetotartrate solu...	20150B	2
Hedinger	Aluminium acetotartrate solu...	20150B	4
Noweda	Aluminium acetotartrate solu...	20150B	1
Hedinger	Aluminium acetotartrate solu...	20150B	1
Hedinger	Aluminium acetotartrate solu...	20150b	1
Stumpf	Aluminium acetotartrate solu...	20151B	1
Kehr	Aluminium acetotartrate solu...	20151B	1
Hedinger	Aluminium acetotartrate solu...	20151B	1
Hedinger	Aluminium acetotartrate solu...	20151B	8
hedinger	Aluminium acetotartrate solu...	20151B	1
Anzag/	Aluminium acetotartrate solu...	20151B	1
HEDINGER	Aluminium acetotartrate solu...	20151B	1
Großhandel	Aluminium acetotartrate solu...	20151B	1
Caelo	Aluminium acetotartrate solu...	20151B	1
Noweda	Aluminium acetotartrate solu...	20151B	1
Hedinger/Noweda	Aluminium acetotartrate solu...	20151B	1
Franz GmbH/Phönix	Aluminium acetotartrate solu...	20151B	1
PHÖNIX 05.06.2013	Aluminium acetotartrate solu...	20151B	1
Dr.Seeger GmbH	Aluminium acetotartrate solu...		1
hedinger	Aluminium acetotartrate solu...	20152B	1
Hedinger/Phönix	Aluminium acetotartrate solu...	20152B	2
Anzag	Aluminium acetotartrate solu...	20152B	2
Hedinger	Aluminium acetotartrate solu...	20152B	6
Sanacorp	Aluminium acetotartrate solu...	20152B	1
Caelo	Aluminium acetotartrate solu...	20152b	1
Anzag 17.04.13 ek: 5,43EUR	Aluminium acetotartrate solu...	20152B	1
Caelo	Aluminium acetotartrate solu...	20152B	2
Hedinger / Gehe	Aluminium acetotartrate solu...	20152B	1
Fischar/Anzag	Aluminium acetotartrate solu...	20152B	1
Hedinger/Gehe	Aluminium acetotartrate solu...	20152B	1
Caelo	Aluminium acetotartrate solu...	20153B	5
Hedinger/Otto	Aluminium acetotartrate solu...	20153B	1
Hedinger	Aluminium acetotartrate solu...	20153B	2
franz gmbh	Aluminium acetotartrate solu...	20153B	1
Spangr/Hedinger	Aluminium acetotartrate solu...	20153B	1
Gehe/Hedinger	Aluminium acetotartrate solu...	20153B	1
Hedinger , Noweda	Aluminium acetotartrate solu...	20153B	1
Hedinger	Aluminium acetotartrate solu...	20153B	10
Hedinger/Phönix	Aluminium acetotartrate solu...	20153B	2
Hedinger/Sanacorp	Aluminium acetotartrate solu...	20153B	2
Gehe	Aluminium acetotartrate solu...	20153B	1
Caelo	Aluminium acetotartrate solu...	20153B	3
hedinger/anzag	Aluminium acetotartrate solu...	20153b	1
Dr.Seeger GmbH	Aluminium acetotartrate solu...	20153B	1
Phönix	Aluminium acetotartrate solu...	20153b	1
Dr. Seeger	Aluminium acetotartrate solu...	20153B	1
Noweda/Hedinger	Aluminium acetotartrate solu...	20153B	1
Hedinger	Aluminium acetotartrate solu...	20154B	1
geile	Aluminium acetotartrate solu...	20154B	1
Caelo	Aluminium acetotartrate solu...	20154B	5
hedinger/gehe	Aluminium acetotartrate solu...	20154B	1
Hedinger/V.d.L.	Aluminium acetotartrate solu...	20154B	1
Hedinger/ Sanacorp	Aluminium acetotartrate solu...	20154B	1
Hedinger, Noweda	Aluminium acetotartrate solu...	20154B	2
Hedinger , Stuttgart	Aluminium acetotartrate solu...	20154B	1
Hedinger/Sanacorp	Aluminium acetotartrate solu...	20154B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
hedinger/Gehe	Aluminium acetotartrate solu. . .	20155B	1
hedinger/phönix	Aluminium acetotartrate solu. . .	20155B	1
Hedinger/GEHE	Aluminium acetotartrate solu. . .	20155B	1
heidinger/noweda	Aluminium acetotartrate solu. . .	20156B	1
Caelo	Aluminium acetotartrate solu. . .	20158B	1
Hedinger	Aluminium acetotartrate solu. . .	2015B	1
Hedinger/Phönix	Aluminium acetotartrate solu. . .	2015B	1
Hedinger	Aluminium acetotartrate solu. . .	2021302	1
Hedinger	Aluminium acetotartrate solu. . .	2153B	1
Hedinger	Aluminium acetotartrate solu. . .	402211793	1
athenstaedt GmbH	Aluminium acetotartrate solu. . .	900107	2
Phönix Gotha	Aluminium acetotartrate solu. . .	900108	2
Athenstaedt	Aluminium acetotartrate solu. . .	900108	3
athenstedt	Aluminium acetotartrate solu. . .	110029	1
Hedinger	Aluminium acetotartrate solu. . .	1111	1
Caelo	Aluminium acetotartrate solu. . .	11315004	1
Caelo	Aluminium acetotartrate solu. . .	11315002	1
Caelo	Aluminium acetotartrate solu. . .	11315007	1
Caelo	Aluminium acetotartrate solu. . .	12116803	3
Caelo	Aluminium acetotartrate solu. . .	12116804	5
Caesar&loretz,Phoenix	Aluminium acetotartrate solu. . .	12116804	1
Caelo	Aluminium acetotartrate solu. . .	12233501	4
Caelo	Aluminium acetotartrate solu. . .	12233502	4
Caelo	Aluminium acetotartrate solu. . .	12233504	1
Caelo	Aluminium acetotartrate solu. . .	12233503	1
Caelo	Aluminium acetotartrate solu. . .	12318701	6
Caelo	Aluminium acetotartrate solu. . .	12318702	1
Caelo	Aluminium acetotartrate solu. . .	12318704	3
Phönix	Aluminium acetotartrate solu. . .	12318703	1
Caelo	Aluminium acetotartrate solu. . .	12318705	2
Phönix	Aluminium acetotartrate solu. . .	12318705	1
Caelo	Aluminium acetotartrate solu. . .	13101101	2
Caelo	Aluminium acetotartrate solu. . .	13101102	5
Phönix	Aluminium acetotartrate solu. . .	13101101	1
Caelo	Aluminium acetotartrate solu. . .	13101104	6
Calo/AHD	Aluminium acetotartrate solu. . .	13101105	1
Caelo	Aluminium acetotartrate solu. . .	13261301	3
Caesar&Loretz GmbH/Anzag	Aluminium acetotartrate solu. . .	13101105	1
Caelo	Aluminium acetotartrate solu. . .	13261306	1
Caelo	Aluminium acetotartrate solu. . .	13261307	1
Caelo	Aluminium acetotartrate solu. . .	17081201	1
Hedinger Phönix	Aluminium acetotartrate solu. . .	201528	1
Hedinger	Aluminium acetotartrate solu. . .	26041301	1
Caelo	Aluminium acetotartrate solu. . .	2713M-01793	1
Hedinger	Aluminium acetotartrate solu. . .		3
Serumwerk Bernburg AG	Water		2
Caelo	Water	15068702	6
Caelo	Water	15287302	3
Caelo	Water	15336501	7
Caelo	Water	15394202	8
Fagron	Water	15F24-T01-023879	9
Fagron	Water	15H25-T14-025670	8
Fagron	Water	15G09-T02-024364	4
Fagron	Water	15K10-T04-027614	5
Caelo	Water	15I10-T08-025964	1
Caelo	Water	15H25-T18-025678	1
Caelo	Water	16127901	3
Fagron	Water	16F13-T11-034264	8
Fagron	Water	16F13-T16-034297	2
Fagron	Water	16F13-T14-034274	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Water	16F14-T07-034326	8
Fagron	Water	16L12-T11-039624	3
Fagron	Water	16L12T16039628	1
Fagron	Water	17E05-T01-044045	2
Fagron	Water	17E05-T05-044063	2
Wepa	Water	LOTAP0025	1
Caelo	Chlorhexidine digluconate so...		1
Fagron	Chlorhexidine digluconate so...	10031503	1
Caelo	Chlorhexidine digluconate so...	10213458	1
Fagron	Chlorhexidine digluconate so...	1081403	1
Caelo	Chlorhexidine digluconate so...	10213441	1
Euro OTC	Chlorhexidine digluconate so...	1206012-02	1
Fagron	Chlorhexidine digluconate so...	11081502	1
Euro OTC	Chlorhexidine digluconate so...	1206023-02	26
EuRho	Chlorhexidine digluconate so...	1206023-02	2
Hedinger/Sanacorp6.5.2014	Chlorhexidine digluconate so...	1234	1
Fragon	Chlorhexidine digluconate so...	12F01-N21	2
Fagron	Chlorhexidine digluconate so...	12F01-N21	40
Euro OTC	Chlorhexidine digluconate so...	1308002	2
Fagron	Chlorhexidine digluconate so...	1308002-001	1
Euro OTC	Chlorhexidine digluconate so...	1308002-01	27
Fagron	Chlorhexidine digluconate so...	12K22-N09	20
Caelo	Chlorhexidine digluconate so...	1308002	1
VDL	Chlorhexidine digluconate so...	130800201	1
Euro OTC	Chlorhexidine digluconate so...	130800201	4
Noweda	Chlorhexidine digluconate so...	1308002-01	2
Gehe	Chlorhexidine digluconate so...	1308002-01	1
Caelo	Chlorhexidine digluconate so...	1308002-01	1
Euro OTC	Chlorhexidine digluconate so...	1308002-02	29
EuRho	Chlorhexidine digluconate so...	1308002-02	1
Caelo	Chlorhexidine digluconate so...	1308002-02	8
Gehe	Chlorhexidine digluconate so...	1308002-02	1
Evonic	Chlorhexidine digluconate so...	1308002-02	2
	Chlorhexidine digluconate so...	1308002-02	1
Caelo	Chlorhexidine digluconate so...	13081502	1
Caelo	Chlorhexidine digluconate so...	13136105	5
Caelo	Chlorhexidine digluconate so...	13136114	10
GEHE 20.05.2014	Chlorhexidine digluconate so...	13136114	1
Berkel AHK	Chlorhexidine digluconate so...	13136115	1
AHD	Chlorhexidine digluconate so...	13136115	1
Caelo	Chlorhexidine digluconate so...	13136115	2
Fagron	Chlorhexidine digluconate so...	13D16-N04	9
Fagron	Chlorhexidine digluconate so...	13A29-N06	14
Fagron	Chlorhexidine digluconate so...	13H19-B01	5
Fagron	Chlorhexidine digluconate so...	13H19-B01-285448	3
Fagron	Chlorhexidine digluconate so...	13H19-B01-285451	6
Fagron	Chlorhexidine digluconate so...	13H19-B01-285449	1
Fagron	Chlorhexidine digluconate so...	13H19-B01-285452	4
Fagron	Chlorhexidine digluconate so...	13H19-B01-285453	9
Fagron	Chlorhexidine digluconate so...	13K16B06288853	1
Fagron	Chlorhexidine digluconate so...	13K26-B06	8
Fagron	Chlorhexidine digluconate so...	13k26-b06	1
Caelo	Chlorhexidine digluconate so...	13K26-B06	1
Fagron	Chlorhexidine digluconate so...	13K26-B06289365	1
Fagron	Chlorhexidine digluconate so...	13K26-B06-288848	2
Caelo	Chlorhexidine digluconate so...	13K26-B06-288849	2
Fagron	Chlorhexidine digluconate so...	13K26-B06-288849	1
PHOENIX 25.07.2014	Chlorhexidine digluconate so...	13k26-b06-288849	1
Fagron	Chlorhexidine digluconate so...	13K26-B06-288853	3
	Chlorhexidine digluconate so...	13K26B06289365	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Chlorhexidine digluconate so...	13K-26-B06-288853	1
Phönix Gotha	Chlorhexidine digluconate so...	13K26-B06-288853	1
Hedinger/Sanacorp6.5.2014	Chlorhexidine digluconate so...	13K26-B06289365	1
Audor Pharma	Chlorhexidine digluconate so...	13K26-B06-289365	2
Fagron	Chlorhexidine digluconate so...	13K26-B06-289365	6
Fagro/ Noweda	Chlorhexidine digluconate so...	13K26-B06-289365	1
Fagron	Chlorhexidine digluconate so...	13K26B06289366	1
Fagron	Chlorhexidine digluconate so...	13K26-B06-289366	1
Caelo	Chlorhexidine digluconate so...	14011004	2
Euro OTC	Chlorhexidine digluconate so...	13KB062888448	1
Caelo	Chlorhexidine digluconate so...	14011009	1
Sanacorp	Chlorhexidine digluconate so...	14011004	1
Caelo	Chlorhexidine digluconate so...	14011013	5
Caelo	Chlorhexidine digluconate so...	14011011	3
sanacorp	Chlorhexidine digluconate so...	14011015	1
Caelo	Chlorhexidine digluconate so...	14011015	5
Lamotte / Gehe	Chlorhexidine digluconate so...	14011017	1
Caelo	Chlorhexidine digluconate so...	14011017	4
Caelo	Chlorhexidine digluconate so...	14011020	1
Caelo	Chlorhexidine digluconate so...	14011018	3
Caelo	Chlorhexidine digluconate so...	1401105	1
Euro OTC	Chlorhexidine digluconate so...	1412018-01	24
Caelo	Chlorhexidine digluconate so...	140111009	1
EuRho	Chlorhexidine digluconate so...	141201801	1
EuroPTC/Phönix	Chlorhexidine digluconate so...	1412018-01	1
Caelo	Chlorhexidine digluconate so...	1412018-01	4
Phönix Gotha	Chlorhexidine digluconate so...	1412018-01	1
VDL; 13.7.15; 3,83EUR	Chlorhexidine digluconate so...	1412018-01	1
Euro OTC	Chlorhexidine digluconate so...	141201801	1
Euro OTC	Chlorhexidine digluconate so...	1412018-02	26
Caelo	Chlorhexidine digluconate so...	141201801	1
Caelo	Chlorhexidine digluconate so...	141201802	2
Fagro/AHD	Chlorhexidine digluconate so...	1412018-02	1
Caelo	Chlorhexidine digluconate so...	1412018-02	2
Anzag; 3.80EUR; 7.10.2015	Chlorhexidine digluconate so...	1412018-02	1
Eu Rho/ Gehe	Chlorhexidine digluconate so...	1412018-02	1
EuRho	Chlorhexidine digluconate so...	1412018-02	1
Bombastus	Chlorhexidine digluconate so...	1412018-02	1
Phönix Gotha	Chlorhexidine digluconate so...	14120189-01	1
Euro OTC	Chlorhexidine digluconate so...	1412018801	1
Caelo	Chlorhexidine digluconate so...	14331301	7
Euro OTC	Chlorhexidine digluconate so...	141218-01	1
Caelo	Chlorhexidine digluconate so...	14331303	2
Noweda18.12.2014	Chlorhexidine digluconate so...	14331303	1
Caelo	Chlorhexidine digluconate so...	14331306	4
Caelo	Chlorhexidine digluconate so...	14331307	8
Phoenix	Chlorhexidine digluconate so...	14331307	1
Cealo	Chlorhexidine digluconate so...	14331307	1
Caelo	Chlorhexidine digluconate so...	14331308	2
Caelo	Chlorhexidine digluconate so...	14331309	3
Caelo	Chlorhexidine digluconate so...	14331312	8
Caelo	Chlorhexidine digluconate so...	14331310	1
sanacorp	Chlorhexidine digluconate so...	14331312	1
Gehe	Chlorhexidine digluconate so...	14331309	1
Sanacorp	Chlorhexidine digluconate so...	14331312	1
Caelo	Chlorhexidine digluconate so...	14331313	9
Fagron	Chlorhexidine digluconate so...	14D25-B04	2
Caelo	Chlorhexidine digluconate so...	14D25-B04	2
Fagron	Chlorhexidine digluconate so...	14D25-B04-294673	4
Fagron	Chlorhexidine digluconate so...	14332309	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Chlorhexidine digluconate so...	14D25-B04-294674	3
Fagron	Chlorhexidine digluconate so...	14D25-B04-294675	1
Fagron	Chlorhexidine digluconate so...	14D25-B04-299289	1
Fagron	Chlorhexidine digluconate so...	14D25-B04-299288	1
Fagron	Chlorhexidine digluconate so...	15/03-B02-315212	1
Fagron	Chlorhexidine digluconate so...	15/03-B02-315213	2
Caelo	Chlorhexidine digluconate so...	15/03-B02-315213	1
Caelo	Chlorhexidine digluconate so...	15/03-B02-315552	1
Euro OTC	Chlorhexidine digluconate so...	15008022-01	1
Fagron	Chlorhexidine digluconate so...	1503B02316428	1
Fagron	Chlorhexidine digluconate so...	1503-B02-316428	1
Euro OTC	Chlorhexidine digluconate so...	1508022-01	24
Caelo	Chlorhexidine digluconate so...	1508022-01	3
Euro OTC	Chlorhexidine digluconate so...	150802201	3
Noweda	Chlorhexidine digluconate so...	1508022-01	1
Apo-ident	Chlorhexidine digluconate so...	1508022-01	1
Euro OTC	Chlorhexidine digluconate so...	1508022-02	10
Caelo	Chlorhexidine digluconate so...	1508022-22	1
Fagron	Chlorhexidine digluconate so...	15103-B02-315213	1
Bombastus	Chlorhexidine digluconate so...	1508022-02	1
Euro OTC	Chlorhexidine digluconate so...	1508022-22	1
Euro OTC	Chlorhexidine digluconate so...	1508072-02	1
Caelo	Chlorhexidine digluconate so...	15109201	11
Caelo	Chlorhexidine digluconate so...	15109202	5
Euro OTC	Chlorhexidine digluconate so...	151103-01	1
Caelo	Chlorhexidine digluconate so...	1511033/02	1
Fagron	Chlorhexidine digluconate so...	1511033-01	2
Euro OTC	Chlorhexidine digluconate so...	1511033-01	44
Gehe	Chlorhexidine digluconate so...	151103301	1
phönix	Chlorhexidine digluconate so...	1511033-01	2
Euro OTC	Chlorhexidine digluconate so...	151103301	1
Caelo	Chlorhexidine digluconate so...	1511033-01	2
Phoenix	Chlorhexidine digluconate so...	1511033-01	1
Caelo	Chlorhexidine digluconate so...	1511033-02	2
Euro OTC	Chlorhexidine digluconate so...	1511033-02	6
Galderma	Chlorhexidine digluconate so...	1511033-02	1
Fagron	Chlorhexidine digluconate so...	15121406	1
Caelo	Chlorhexidine digluconate so...	15131502	5
Caelo	Chlorhexidine digluconate so...	15131503	11
Geilenkirchen	Chlorhexidine digluconate so...	15131503	1
anzag	Chlorhexidine digluconate so...	15131504	1
Caelo	Chlorhexidine digluconate so...	15131504	10
Sanacorp 17.07.2015	Chlorhexidine digluconate so...	15131504	1
Caelo	Chlorhexidine digluconate so...	15131505	14
Fagron	Chlorhexidine digluconate so...	15131505	2
Caelo	Chlorhexidine digluconate so...	15131508	2
Caelo	Chlorhexidine digluconate so...	15131509	21
Euro OTC	Chlorhexidine digluconate so...	15131505	1
Sanacorp	Chlorhexidine digluconate so...	15131509	1
Caelo	Chlorhexidine digluconate so...	15131510	12
Noweda	Chlorhexidine digluconate so...	15131510	1
Caelo	Chlorhexidine digluconate so...	15131512	7
Caelo	Chlorhexidine digluconate so...	15131604	1
Caelo	Chlorhexidine digluconate so...	151361503	1
Caelo	Chlorhexidine digluconate so...	15292504	16
Frau Brinn-Holtz	Chlorhexidine digluconate so...	15292507	1
Caleo/Sanacorp	Chlorhexidine digluconate so...	15292507	1
Caelo	Chlorhexidine digluconate so...	15292507	5
Caelo	Chlorhexidine digluconate so...	15292508	11
Sanacorp	Chlorhexidine digluconate so...	15292508	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Chlorhexidine digluconate so...	15292508	1
Caesar & Loretz GmbH/ A...	Chlorhexidine digluconate so...	15292508	1
Cealo/Noweda	Chlorhexidine digluconate so...	15292508	2
Caelo	Chlorhexidine digluconate so...	15292509	6
Caelo	Chlorhexidine digluconate so...	15292511	2
Caelo	Chlorhexidine digluconate so...	15292512	4
Caelo	Chlorhexidine digluconate so...	1529504	1
Caelo	Chlorhexidine digluconate so...	1529509	1
Euro OTC	Chlorhexidine digluconate so...	1529508	1
Caelo	Chlorhexidine digluconate so...	1529511	1
Caelo	Chlorhexidine digluconate so...	1592509	1
Fagron	Chlorhexidine digluconate so...	15C03B02	2
Fagron	Chlorhexidine digluconate so...	15C03-B02-306308	2
Fragon/ Phönix	Chlorhexidine digluconate so...	15C03-B02-306308	1
Schütz Dishman Biotech ...	Chlorhexidine digluconate so...	15E13-B07	1
Fagron	Chlorhexidine digluconate so...	15E13-B07-309826	5
Fagron	Chlorhexidine digluconate so...	15E13B07309826	1
Fagron	Chlorhexidine digluconate so...	15F30-B05	3
Fagron	Chlorhexidine digluconate so...	15F30B05	1
Fagron	Chlorhexidine digluconate so...	15F30-B05-309827	1
Fagron	Chlorhexidine digluconate so...	15F30-B05-309829	2
Fagron	Chlorhexidine digluconate so...	15F30-B05-311700	4
Caelo	Chlorhexidine digluconate so...	15f30-b05-309829	1
Caelo	Chlorhexidine digluconate so...	15F30-B05-311700	1
Fagron	Chlorhexidine digluconate so...	15I03-B01-314176	1
Fagron	Chlorhexidine digluconate so...	15I03-B02	1
Fagron	Chlorhexidine digluconate so...	15I03B02	1
Fagron	Chlorhexidine digluconate so...	15I03-B02-0006	1
Fagron	Chlorhexidine digluconate so...	15I03-B02-313663	4
Schütz Dishman Biotech ...	Chlorhexidine digluconate so...	15I03-B02	1
Fagron	Chlorhexidine digluconate so...	15I03B02314176	1
Caelo	Chlorhexidine digluconate so...	15I03-B02-315213	4
Fagron	Chlorhexidine digluconate so...	15I03-B02-314177	1
Fagron	Chlorhexidine digluconate so...	15I03-B02-315213	1
Fagron	Chlorhexidine digluconate so...	15I03-B02-315552	4
Fagron	Chlorhexidine digluconate so...	15I03B02317319	2
Fagron	Chlorhexidine digluconate so...	15L03-B02-315552	1
Fagron	Chlorhexidine digluconate so...	15J03-B02	1
Caelo	Chlorhexidine digluconate so...	15I03-b02-317319	1
Sonnen-Apotheke Wefensl...	Chlorhexidine digluconate so...	16002501	1
Fagron	Chlorhexidine digluconate so...	15I03-B02-313663	1
Caelo	Chlorhexidine digluconate so...	16002501	10
Sanacorp	Chlorhexidine digluconate so...	16002501	2
Caelo	Chlorhexidine digluconate so...	16002502	2
Caelo	Chlorhexidine digluconate so...	16002505	10
Caelo	Chlorhexidine digluconate so...	16002506	4
Ichthyol/Phönix	Chlorhexidine digluconate so...	16002506	2
Caesar & Loretz	Chlorhexidine digluconate so...	16002510	1
Caelo	Chlorhexidine digluconate so...	16002510	5
Fagron	Chlorhexidine digluconate so...	16002510	2
hochhaus-Apotheke	Chlorhexidine digluconate so...	16002510	1
Euro OTC	Chlorhexidine digluconate so...	1610021-01	10
Evonik Technochem. GmbH	Chlorhexidine digluconate so...	1610021-01	1
Phönix	Chlorhexidine digluconate so...	1610021-01	1
Euro/Phönix	Chlorhexidine digluconate so...	161002101	2
Euro OTC	Chlorhexidine digluconate so...	1610022-01	5
Caelo	Chlorhexidine digluconate so...	16308801	23
GEHE	Chlorhexidine digluconate so...	1610022-01	1
Gehe	Chlorhexidine digluconate so...	16308801	1
Caelo	Chlorhexidine digluconate so...	16308802	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chlorhexidine digluconate so...	16308811	13
Caelo	Chlorhexidine digluconate so...	16308812	3
Phönix	Chlorhexidine digluconate so...	16308812	1
Fagron	Chlorhexidine digluconate so...	16B25-B98	1
Fagron	Chlorhexidine digluconate so...	16B25-B98-318800	4
Caelo	Chlorhexidine digluconate so...	16309902	1
Fagron	Chlorhexidine digluconate so...	16B25-B98-318799	1
Fagron	Chlorhexidine digluconate so...	16e09-b03	1
Fagron	Chlorhexidine digluconate so...	16E09-B03	2
Caelo	Chlorhexidine digluconate so...	16E09-B03	1
Fagron	Chlorhexidine digluconate so...	16E09-B03-320995	5
Fagron	Chlorhexidine digluconate so...	16E09-B03-320996	2
Fagron	Chlorhexidine digluconate so...	16E09-B03-322382	2
Fagron	Chlorhexidine digluconate so...	16E09-B03-324024	1
Fagron	Chlorhexidine digluconate so...	16H-23-B03-329123	1
Fagron	Chlorhexidine digluconate so...	16E09-B03-322464	2
Caelo	Chlorhexidine digluconate so...	16E09-B03-324024	1
Fagron	Chlorhexidine digluconate so...	16H23-B03-329123	1
Fagron	Chlorhexidine digluconate so...	16E09-B03-323873	1
Phönix	Chlorhexidine digluconate so...	17001037-02	1
Sanacorp EK:5,70 EUR...	Chlorhexidine digluconate so...	1701037	1
Euro OTC	Chlorhexidine digluconate so...	1701037	2
Euro OTC	Chlorhexidine digluconate so...	1701037-01	5
Finteler Apotheke	Chlorhexidine digluconate so...	1701037-01	1
Phoenix WE: 01.09.17 E...	Chlorhexidine digluconate so...	1701037-02	3
Euro OTC	Chlorhexidine digluconate so...	1701037-02	6
Fagron	Chlorhexidine digluconate so...	1701037-01	1
Caelo	Chlorhexidine digluconate so...	1701037-02	1
Euro OTC	Chlorhexidine digluconate so...	1708006-01	4
Caelo	Chlorhexidine digluconate so...	170103702	1
Euro OTC	Chlorhexidine digluconate so...	1708006-02	3
Caelo	Chlorhexidine digluconate so...	17210101	1
Caelo	Chlorhexidine digluconate so...	171027CDG/17210105	1
Caelo	Chlorhexidine digluconate so...	17210105	7
Caelo	Chlorhexidine digluconate so...	17210106	3
Noweda	Chlorhexidine digluconate so...	17210105	1
Caelo	Chlorhexidine digluconate so...	17210110	3
GEHE	Chlorhexidine digluconate so...	17210110	1
Noweda	Chlorhexidine digluconate so...	17210112	2
Allianz Healthcare	Chlorhexidine digluconate so...	1721110	1
Caelo	Chlorhexidine digluconate so...	17230801	11
EuRho/Noweda	Chlorhexidine digluconate so...	17230801	1
Fagron	Chlorhexidine digluconate so...	17E09-B01-335136	2
Fagron	Chlorhexidine digluconate so...	17E09-B01-335140	1
Caelo	Chlorhexidine digluconate so...	17E09-B01-335140	1
Fagron	Chlorhexidine digluconate so...	17e09-b01-335140	1
Fagron	Chlorhexidine digluconate so...	18A09-B05-344604	1
Fagron	Chlorhexidine digluconate so...	26091402	1
Fagron	Chlorhexidine digluconate so...	2614Q-01262	1
Fagron	Chlorhexidine digluconate so...	18A09-B05-344605	1
Evonik	Chlorhexidine digluconate so...	3896220404	1
Euro OCT	Chlorhexidine digluconate so...	3-8962-3-01-24	1
Euro OTC	Chlorhexidine digluconate so...	3896250214	1
Euro OTC	Chlorhexidine digluconate so...	3-8962-6-03-21	1
Euro OTC	Chlorhexidine digluconate so...	3896270218	2
Euro OTC	Chlorhexidine digluconate so...	3896260321	1
Euro OTC	Chlorhexidine digluconate so...	3-8962-7-02-39	1
Fagron	Chlorhexidine digluconate so...	BK13-CG-018	2
Fagron	Chlorhexidine digluconate so...	5x15I03-B02-3152120	1
Fagron	Chlorhexidine digluconate so...	BK13-CG-285447	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Chlorhexidine digluconate so...	BK14-CG-003	2
GEHE 08.12.2014	Chlorhexidine digluconate so...	BK14-CG-003	1
Fagron	Chlorhexidine digluconate so...	BK14-CG-023	1
Fagron	Chlorhexidine digluconate so...	BK15-CG-001	3
Fagron	Chlorhexidine digluconate so...	BK15CG001	1
Fagron	Chlorhexidine digluconate so...	BK16-CG-005	1
Euro OTC	Chlorhexidine digluconate so...	s-8962-6-03-21	1
Fagron	Chlorhexidine digluconate so...	103311262	1
Fagron	Chlorhexidine digluconate so...	1112013D	1
Fagron	Chlorhexidine digluconate so...	4111303	1
Fagron	Chlorhexidine digluconate so...	6111210	1
Fagron	Chlorhexidine digluconate so...	204211262	1
Caelo	Chlorhexidine digluconate so...	10200303	1
Fagron	Chlorhexidine digluconate so...	BK16-CG-011	1
Fagron	Chlorhexidine digluconate so...	9081302	1
Euro OTC	Chlorhexidine digluconate so...	106023-02	1
Euro OTC	Chlorhexidine digluconate so...	1105012-01	2
EuRho	Chlorhexidine digluconate so...	1105012-01	1
Celo	Chlorhexidine digluconate so...	110816-7	1
Gehe	Chlorhexidine digluconate so...	1110012-01	1
Euro OTC	Chlorhexidine digluconate so...	1110012-01	1
Caelo	Chlorhexidine digluconate so...	11135114	2
Phönix, 17.10.2012	Chlorhexidine digluconate so...	1110012-01	1
Caelo	Chlorhexidine digluconate so...	11135107	1
Gehe	Chlorhexidine digluconate so...	11135114	2
Noweda	Chlorhexidine digluconate so...	11353004	1
Caelo	Chlorhexidine digluconate so...	11353001	3
Caelo	Chlorhexidine digluconate so...	11353004	4
Phönix	Chlorhexidine digluconate so...	11353004	1
Caelo	Chlorhexidine digluconate so...	11353009	2
Caelo	Chlorhexidine digluconate so...	11353012	3
ANZAG	Chlorhexidine digluconate so...	11353012	1
Anzag 21.12.12	Chlorhexidine digluconate so...	11353012	2
Caelo	Chlorhexidine digluconate so...	11353013	5
Noweda	Chlorhexidine digluconate so...	11353013	1
Fagron	Chlorhexidine digluconate so...	11K09-N05	4
Fagron	Chlorhexidine digluconate so...	11L05-N08	5
Caelo	Chlorhexidine digluconate so...	11L05-N08	1
Fagron	Chlorhexidine digluconate so...	120201001	1
Euro OTC	Chlorhexidine digluconate so...	1202010-01	10
Phoenix 17.07.2012	Chlorhexidine digluconate so...	1202010-01	2
Fagron	Chlorhexidine digluconate so...	11L05-N080069/262	1
Noweda	Chlorhexidine digluconate so...	1202010-01	3
Phönix/EuRho	Chlorhexidine digluconate so...	1206023-01	2
Fagron	Chlorhexidine digluconate so...	120355-ooo4	1
Alliance Healthcare	Chlorhexidine digluconate so...	1206023	1
Sanacorp	Chlorhexidine digluconate so...	1206023-01	1
Euro OTC	Chlorhexidine digluconate so...	120602301	3
Euro OTC	Chlorhexidine digluconate so...	1206023-01	15
Phönix	Chlorhexidine digluconate so...	1206023-01	1
Caelo	Chlorhexidine digluconate so...	1206023-01	1
Anzag	Chlorhexidine digluconate so...	1206023-01	1
10 05.06.2013	Chlorhexidine digluconate so...	1206023-01	1
Gehe	Chlorhexidine digluconate so...	1206023-01	1
Gehe	Chlorhexidine digluconate so...	1206023-02	7
Caelo	Chlorhexidine digluconate so...	1206023-02	3
Phönix	Chlorhexidine digluconate so...	1206023-02	2
Phönix 28.01.2014	Chlorhexidine digluconate so...	1206023-02	2
Caelo	Chlorhexidine digluconate so...	12114001	5
Noweda	Chlorhexidine digluconate so...	12114001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Chlorhexidine digluconate so...	12114001	1
Caelo	Chlorhexidine digluconate so...	12114007	1
Fagron	Chlorhexidine digluconate so...	12114007	1
Caelo	Chlorhexidine digluconate so...	12114008	5
Sanacorp	Chlorhexidine digluconate so...	12114001	1
Caelo	Chlorhexidine digluconate so...	12114012	12
Gehe	Chlorhexidine digluconate so...	12114012	1
Noweda	Chlorhexidine digluconate so...	12114012	1
AHD	Chlorhexidine digluconate so...	12114012	1
Caelo	Chlorhexidine digluconate so...	12114013	2
Caelo	Chlorhexidine digluconate so...	12114014	5
Gehe	Chlorhexidine digluconate so...	12114014	1
Caelo	Chlorhexidine digluconate so...	12114019	3
Caelo	Chlorhexidine digluconate so...	12114015	2
Fagron	Chlorhexidine digluconate so...	1213E-01262	1
Fagron	Chlorhexidine digluconate so...	12A02-N11	6
Fagron	Chlorhexidine digluconate so...	12a02-n11	1
Phönix Gotha	Chlorhexidine digluconate so...	12A02-N11	1
Fagron	Chlorhexidine digluconate so...	12F01-N021	1
Fargon	Chlorhexidine digluconate so...	12F01-N21	1
Großhandel	Chlorhexidine digluconate so...	12F01-N21	1
Fargon/ Noweda	Chlorhexidine digluconate so...	12f01n21	1
Fagron	Chlorhexidine digluconate so...	12f01-n21	1
Caelo	Chlorhexidine digluconate so...	12H14-N016	1
Fagron	Chlorhexidine digluconate so...	12H14-N06	22
Caelo	Chlorhexidine digluconate so...	12H14-N06	2
Fagron	Chlorhexidine digluconate so...	12H-14-N06	1
Phönix	Chlorhexidine digluconate so...	12H14-N06	1
Phönix Gotha	Chlorhexidine digluconate so...	12H14-N06	1
Gehe	Chlorhexidine digluconate so...	12h14-n06	1
Fagron	Chlorhexidine digluconate so...	12H14N06	1
ANZAG 27.12.2012	Chlorhexidine digluconate so...	12H14No6	1
Fagron	Chlorhexidine digluconate so...	12H14-No6	4
Fagron	Chlorhexidine digluconate so...	12H14-No06	1
Sanacorp	Chlorhexidine digluconate so...	12K22-N09	1
Fagron	Chlorhexidine digluconate so...	12K22-0005	1
Caelo	Chlorhexidine digluconate so...	12k22-n09	1
Euro OTC	Chlorhexidine digluconate so...	12K22-N09	1
Fagron	Chlorhexidine digluconate so...	12K22-No9	1
Euro OTC	Chlorhexidine digluconate so...	130002-01	1
Caelo	Chlorhexidine digluconate so...	13136015	1
Fagron	Chlorhexidine digluconate so...	13081204	1
Caelo	Chlorhexidine digluconate so...	13136101	4
Caelo	Chlorhexidine digluconate so...	13136103	16
Sanacorp	Chlorhexidine digluconate so...	13136103	1
Caelo	Chlorhexidine digluconate so...	13136112	4
Gehe	Chlorhexidine digluconate so...	13136105	1
Gehe	Chlorhexidine digluconate so...	13136109	1
Cael/Nowda	Chlorhexidine digluconate so...	13136115	1
Sanacorp	Chlorhexidine digluconate so...	13136115	2
Caelo	Chlorhexidine digluconate so...	1316103	1
Caelo	Chlorhexidine digluconate so...	13A29-N06	1
Phönix Gotha	Chlorhexidine digluconate so...	13A29-N06	1
Fagron	Chlorhexidine digluconate so...	13A29-N060155/228	1
Fagron	Chlorhexidine digluconate so...	13D16-N040090/192	1
Klenk	Chlorhexidine digluconate so...	13D16-N04	1
Fagron	Chlorhexidine digluconate so...	13d16-N04	2
Fagron	Chlorhexidine digluconate so...	13H19B01	1
Euro OTC	Chlorhexidine digluconate so...	13D16-N04	1
Fagron	Chlorhexidine digluconate so...	13H19-801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Großhandel	Chlorhexidine digluconate so...	13H19-B01	1
Fagron	Chlorhexidine digluconate so...	13h19b01285447	1
Fagron	Chlorhexidine digluconate so...	13H19-B01-285447	2
Caelo	Chlorhexidine digluconate so...	13H19-B01-285451	1
Fagron	Chlorhexidine digluconate so...	13H19B01285451	1
Phoenix	Chlorhexidine digluconate so...	13H19-B01-285452	1
Fagron	Chlorhexidine digluconate so...	13H19-B01-CG-011	1
Caelo	Chlorhexidine digluconate so...	13k26-b06-28884	1
Caelo	Chlorhexidine digluconate so...	20031301	1
Fagron	Chlorhexidine digluconate so...	20200112-1	1
Fagron	Chlorhexidine digluconate so...	17011310	1
Fagron	Chlorhexidine digluconate so...	15103B02-315212	1
Fagron	Chlorhexidine digluconate so...	23121307	1
Fagron	Chlorhexidine digluconate so...	25051304	1
Fagron	Chlorhexidine digluconate so...	2512M01262	1
Caelo	Chlorhexidine digluconate so...	26A080415	1
Fagron	Chlorhexidine digluconate so...	2612M01262	1
Euro OTC	Chlorhexidine digluconate so...	282215	1
Caelo	Chlorhexidine digluconate so...	32/AD1	1
Euro OTC	Chlorhexidine digluconate so...	3896220404	1
Fagron	Chlorhexidine digluconate so...	4513U1262	1
Caelo	Chlorhexidine digluconate so...	BK13-CG-011	1
Fagron	Chlorhexidine digluconate so...	BK13-CG-011	3
Fiebig	Chlorhexidine digluconate so...	1110012-01	1
Caelo	Chlorhexidine digluconate so...	11L05-N080040	1
Caelo	Chlorhexidine digluconate so...	1202010-01	1
AHD	Chlorhexidine digluconate so...	1206023-02	3
Fagron	Chlorhexidine digluconate so...	12A02-11	1
Phönix	Chlorhexidine digluconate so...	13K26-B06-28965	1
Hedinger / Gehe	Aluminium acetotartrate solu...	20155B	1
Euro OTC	Aluminium acetotartrate solu...	20157B	1
hedinger	Aluminium acetotartrate solu...	20161B	1
heidinger/gehe	Aluminium acetotartrate solu...	201618	1
Caelo	Aluminium acetotartrate solu...	110024	1
Caelo	Aluminium acetotartrate solu...	14252011	1
Caelo	Aluminium acetotartrate solu...	13378702	1
Dr. Seeger GmbH/Noweda	Aluminium acetotartrate solu...	20146B	1
Hedinger/AHD	Aluminium acetotartrate solu...	20146B	1
Caelo	Aluminium acetotartrate solu...	19052014A	1
Hedinger/ Sanacorp	Aluminium acetotartrate solu...	20150B	1
Phönix	Aluminium acetotartrate solu...	20150B	1
hedinger/anzag	Aluminium acetotartrate solu...	20152B	1
Anzag/ Hedinger	Aluminium acetotartrate solu...	20153B	1
Franz GmbH	Aluminium acetotartrate solu...	20153B	1
Hedinger / Gehe	Aluminium acetotartrate solu...	20154B	1
Hedinger/Phoenix	Aluminium acetotartrate solu...	20154B	1
franz GmbH	Aluminium acetotartrate solu...	20155B	1
Caelo	Aluminium acetotartrate solu...	4071201	1
Caelo	Aluminium acetotartrate solu...	10114803	1
Hedinger	Aluminium acetotartrate solu...	1002211793	1
Phönix Gotha	Aluminium acetotartrate solu...	11315007	2
Caelo	Aluminium acetotartrate solu...	12116802	1
Phönix Gotha	Aluminium acetotartrate solu...	12116803	1
Alliance Healthcare	Aluminium acetotartrate solu...	12318705	2
Caelo	Aluminium acetotartrate solu...	13101105	1
Hedinger / Gehe	Aluminium acetotartrate solu...	3958/V09/15	1
Hedinger/Phönix	Aluminium acetotartrate solu...	20155B	1
Hedinger/AHD	Aluminium acetotartrate solu...	20152B	2
Hedinger/Holdermann	Aluminium acetotartrate solu...	20150B	1
Franz GmbH / Phoenix	Aluminium acetotartrate solu...	201538	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hedinger/AHD	Aluminium acetotartrate solu...	201153B	1
Hedinger Phönix	Aluminium acetotartrate solu...	20153B	1
Caesar & Lorenz	Aluminium acetotartrate solu...	12318703	1
Industries AG / Phoenix	Chlorhexidine digluconate so...	1308002-02	1
Fagron	Chlorhexidine digluconate so...	15/03-B02-313663	1
Euro OTC	Chlorhexidine digluconate so...	150822-02	1
apo-ident.	Chlorhexidine digluconate so...	1511033-01	6
apo-ident.	Chlorhexidine digluconate so...	1610021-01	1
Fagron	Chlorhexidine digluconate so...	16b25-b98-318800	1
Euro OTC	Chlorhexidine digluconate so...	17080006-02	1
A.H.D.	Chlorhexidine digluconate so...	1206023-01	3
Caelo	Chlorhexidine digluconate so...	12122222	1
Caelo	Chlorhexidine digluconate so...	1206023-02-335731907	1
Wepa/ wepa	Water		1
Fagron	Water	13K13-T08-008092	1
Caelo	Water	14170717	1
Fagron	Water	14E30-T11-012852	1
Fagron	Water	14F02-T03-012891	1
Caelo	Water	14321201	1
Caelo	Water	14364201	1
Fagron	Water	14H21-T05-014956	1
Fagron	Water	14K05-T02-016775	2
Phönix Gotha	Water	14K05-T02-016775	2
Fagron	Water	14K05-T03-016779	1
Fagron	Water	14J16-T24-016261	1
Fresenius	Water	14KB1003	1
Fagron	Water	14L03-T03-017773	1
Fagron	Water	14L03-T07-017788	3
Fagron	Water	14L03-T09-017794	2
Fagron	Water	14L03-T11-017800	1
Rabenapotheke	Water	14L03-T11-017800	1
Fagron	Water	14L03T09017794	1
Fagron	Water	15/10-T08-025964	1
Fagron	Water	14L03-T13-017804	1
Caelo	Water	15013201	1
Phönix Gotha	Water	15019101	3
Fagron	Water	15/10-T09-025968	1
Caelo	Water	15025501	1
Caelo	Water	15068701	1
Phönix Gotha	Water	15142501	1
Caelo	Water	15220601	3
Caelo	Water	15220602	6
Fagron	Water	1509-T02-024364	1
Phönix Gotha	Water	15235201	2
Caelo	Water	15235301	2
Caelo	Water	1524101	1
Phönix Gotha	Water	15253901	1
Caelo	Water	15253901	3
Caelo	Water	15254101	3
Caelo	Water	15287401	2
Caelo	Water	15287501	2
Caelo	Water	15317401	1
Caelo	Water	153712002	1
Caelo	Water	15355301	1
Caelo	Water	15371202	4
Caelo	Water	15336502	1
Caelo	Water	15355401	1
Caelo	Water	15394201	2
Caesar & Loretz GmbH Hi...	Water	15394202	1
Caelo	Water	15394501	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Water	15427201	5
Caelo	Water	15406801	2
Caelo	Water	1542701	1
Fagron	Water	15D13-T07-021465	2
Fagron	Water	15D13-T09-021469	6
Klenk	Water	15D14-T-021499	1
Fagron	Water	15D14-T05	2
Fagron	Water	15D14-T05-021481	4
Fagron	Water	15D14-T06-021484	4
Wasser, gereinigt	Water	15D14T09021490	1
Fagron	Water	15D14-T09-021490	3
Fagron	Water	15D14-T10-021493	4
Fagron	Water	15D14-T13-021499	2
Fagron	Water	15D14-T14-021502	6
Fargon / Gehe	Water	15D14-T13-021499	1
Fagron	Water	15F04-T03-023259	1
Fagron	Water	15F04-T03-023529	1
Fagron	Water	15D14-To5-021481	1
Noweda	Water	15F04-T03-023259	1
Caelo	Water	15F04-T03-23259	1
Caelo	Water	15F04-T05-023263	1
Fagron	Water	15F04-T05-023263	3
Fargon/Kehr	Water	15F04-T05-023263	1
Fagron	Water	15F-24-023883	1
Fagron	Water	15F24-T03-023883	3
Fagron	Water	15F24T03-023883	1
Caelo	Water	15F24-T03-023883	2
Fagron	Water	15G09-T04-024368	4
Fagron	Water	15h0-T09-025968	1
Caelo	Water	15G09-T04-024368	1
Fagron	Water	15H25-T16-025674	4
Fagron	Water	15H25-T18-025678	4
Fagron	Water	15H25-T21-025683	5
Caelo	Water	15H25-T21-025683	2
Fagron	Water	15h25-t21-025683	1
Fagron	Water	15I10-T08-025964	1
Fagron	Water	15k10-t04-027614	1
Ebert 04.04.2014	Water	15K10-T04-027614	1
Noweda	Water	15K10-T05-027617	1
Caelo	Water	15K10-T06-027620	1
Fargon GmbH & Co.KG	Water	15K10-T07-027623	2
Fagron	Water	15K10-T06-027620	3
Fargon GmbH & Co.KG	Water	15K10-To4-027614	1
Fragon/Phoenix	Water	15L14-T12-028660	1
Fagron	Water	15L14-T12-028660	3
Henry Lamotte	Water	15L14-T12-028660	1
Phoenix	Water	15L14-T14	1
Fagron	Water	15L14-T14-028664	6
Fagron	Water	15L15-T01-028672	6
Fagron	Water	15L15-T01-0208672	1
Fagron	Water	15L14T14028664	2
Fragron	Water	15L15-T01-028672	1
Fargon/Noweda	Water	15L15-T028676	1
Fagron	Water	15L15-T03-028676	2
Fagron	Water	15L15-T03-0286761042	1
Fagron	Water	15L15-T05-028681	2
Caelo	Water	15L15-T05-028681	1
Fagron	Water	15L15-T07-028685	5
Phoenix	Water	15L15-T07	1
Euro OTC	Water	15L15-T07-028685	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Water	15L-15T07028685	1
Fagron	Water	15L15T07-028685	1
Caelo	Water	16006001	3
Caelo	Water	16002901	1
Caelo	Water	16002601	1
Fragon/Phoenix	Water	15L4-T14-028664	1
Alte Apotheke	Water	16022202	1
Caelo	Water	16065001	1
Caelo	Water	16065101	1
Caelo	Water	16065601	1
Caelo	Water	16093301	3
Caesar & Loretz GmbH Hi...	Water	16127902	1
Caelo	Water	16151201	1
Caelo	Water	16127902	1
Phönix	Water	16127901	1
Caelo	Water	1615201	1
Caelo	Water	16178101	1
Caelo	Water	16211901	3
Caelo	Water	16215202	1
Caelo	Water	16223601	3
Caelo	Water	16218601	1
Alte Apotheke	Water	16223601	1
Caelo	Water	16223801	10
Caelo	Water	16255201	2
Caelo	Water	16223901	1
Caelo	Water	16255101	1
Fagron	Water	16286801	2
Caelo	Water	16286801	4
Caelo	Water	16313201	3
Caelo	Water	16307501	1
Caelo	Water	16334301	1
Caelo	Water	16340201	2
Caelo	Water	16313202	1
Euro OTC	Water	16307601	1
Caelo	Water	16353303	4
Fargon GmbH & Co.KG	Water	16B10-T04-030213	3
Fagron	Water	16B10-T02-030206	1
Fagron	Water	16C22-T09-031414	4
Fagron	Water	16B10-T05-03026	1
Fagron	Water	16C23-T06-031423	4
Caelo	Water	16C23-/08031427	1
Fagron	Water	16C23-T06-31423	1
Fagron	Water	16C23-T08-031427	2
WE:Sanacorp 12.09.2016 ...	Water	16C23-T08-031427	1
Fagron	Water	16C23-T14-031432	2
Caelo	Water	16C23-T16-031436	2
Fagron	Water	16C23-T16-031436	2
Fargon GmbH & Co.KG	Water	16F13-T16-034297	1
Fagron	Water	16F13-T18-034301	5
Fagron	Water	16F13T18034301	2
Euro OTC	Water	16F13-T18-034301	1
Fagron	Water	16F14T05034321	1
Fagron	Water	16F14-t03-034317	1
Caelo	Water	16F14-T03-034317	1
Fagron	Water	16F14-T03-034317	4
Fagron	Water	16F14T05-034321	1
Fagron	Water	16F14-T05-034321	4
Fagron	Water	16f14-t07-034326	1
Caelo	Water	16F14-T07-034326	1
Fargon GmbH & Co.KG	Water	16F14-T07-034326	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Water	16H02T01035743	1
Caelo	Water	16K25-C20-17409	1
Phoenix WE28.04.17 -7,9...	Water	16H02-T01-035743	1
Fagron	Water	16F14-T09-034330	2
Fagron	Water	16K25C20-17409	1
Caelo	Water	16L12T16039628	1
Caelo	Water	16L12-T11-039624	1
Fargon	Water	16L12T11039624	1
Fagron	Water	16L12-T16-039628	3
Fagron	Water	16L12-T25-039642	3
Fagron	Water	16L12-T27-039646	2
Fagron	Water	16L12-t29-039657	1
Fagron	Water	16L12-T29-039657	1
Fagron	Water	16L12T29039657	1
Caelo	Water	17002001	4
Fagron	Water	16L1-T27-039646	1
Caelo	Water	17002101	4
Caelo	Water	17032401	2
Caelo	Water	17035201	1
Caelo	Water	17011603	1
Caelo	Water	17034901	1
Caelo	Water	17098701	4
Caelo	Water	17089301	1
Caelo	Water	17098901	5
apomix/Noweda	Water	17114501	1
Caelo	Water	17114501	2
Caelo	Water	17114701	1
Fagron	Water	17114801	1
Caelo	Water	17114801	1
Caelo	Water	17145002	1
Caelo	Water	17145001	1
Klenk	Water	17173001	2
Caelo	Water	17173001	2
Caelo	Water	17179001	3
Caelo	Water	17177301	3
Caelo	Water	17173101	1
Caelo	Water	17179101	3
Caelo	Water	17219101	2
Caelo	Water	17233001	2
Caelo	Water	17233101	3
Caelo	Water	1723311	1
Caelo	Water	17282701	4
Caelo	Water	17295701	3
Caelo	Water	17320001	1
Gehe	Water	17319901	2
Fargon GmbH & Co.KG	Water	17B10-T20-041571	1
Fargon / Gehe	Water	17B10-T27-041575	1
Caelo	Water	17B10-T20-041571	1
Fargon/Kehr	Water	17B10-T20-041571	1
Fagron	Water	17B10-T27-041575	1
GEHE	Water	17B10-T27-041575	2
Fagron	Water	17B13-T02-041604	3
Fagron	Water	17B13-T04-041608	2
Fagron	Water	17C06-T09-042290	1
Sanacorp W.E:16.03.18 ...	Water	17C31-T05-043064	2
Fargon GmbH & Co.KG	Water	17E05-T03-044051	1
Fagron	Water	17E05-T03-044051	1
Fagron	Water	17E05-T13-044080	3
Fagron	Water	17E05-T11-044076	1
Fagron	Water	17F23-T06-045760	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Water	17H24-T09-047709	1
Fagron	Water	17h24-709-047709	1
Fagron	Water	17H14-T09-047709	1
Caelo	Water	17H24-T09-047709	1
Fagron	Water	17H24-T11-047713	2
GEHE	Water	17H24-T13-047717	1
Fargon/Kehr	Water	17H24-T09-047709	1
Fagron	Water	17H24-T15-047721	1
Gehe	Water	17H24-T15-047721	1
GEHE	Water	17H27-T10-050547	1
Fagron	Water	17H24-Tn-047713	1
Fagron	Water	17K01-T01-048440	2
Fagron	Water	17K27-T09-050544	2
Fagron	Water	17L20-T02-051255	4
Caelo	Water	18065601	1
Caelo	Water	18061901	3
Caelo	Water	18093501	1
Fagron	Water	18C13-T14-053757	1
Fagron	Water	219167	1
Caelo	Water	23081606	1
Fargon	Water	303429	1
Hedinger/ Sanacorp	Water	AP0022	1
Caelo	Water	AP0026	1
Wepa	Water	AP0024	1
wepa/Sanacorp	Water	AP0027	1
Fagron	Water	Ap0048	1
WEPA	Water	LOTAP0012	2
Wepa	Water	ap0069	1
Caelo	Water		1
Wepa	Water	LOTAP0015	1
Rondell Apotheke	Water	LOTAP0016	1
Wepa	Water	LOTAP0017	1
Wepa	Water	LOTAP0020	1
WEPA	Water	LOTAP0025	2
Caelo	Water	LOTAP0025	1
wepa/Sanaco	Water	LOTAP0028	1
wepa	Water	LOTAP0033	2
WEPA/Noweda	Water	LOTAP0033	1
Caelo	Water	LOTAP0029	1
Wepa	Water	LOTAP0035	1
WEPA	Water	LOTAP0036	1
wepa	Water	LOTAP0038	1
Wepa/ Noweda	Water	LOTAP0038	1
WEPA	Water	LOTAP0038	2
wepa	Water	LOTAP0040	1
Wepa	Water	LOTAP0040	2
Sanacorp /Wepa	Water	LOTAp0042	1
WEPA	Water	LOTAP0039	1
Caelo	Water	LOTAP0042	1
wepa	Water	LOTAP0045	1
Wepa/Sanacorp	Water	LOTAP0047	1
WEPA Apothekenbedarf	Water	LOTAP0042	1
Wepa / Wepa	Water	LOTAP0045	1
Wepa/ wepa	Water	LOTAP0048	1
Wepa / Wepa	Water	LOTAP0048	2
Fagron	Water	LOTAP0048	1
Caelo	Water	LOTAP0050	1
Wepa/ Wepa	Water	LOTAP0051	1
Wepa/ Phoenix	Water	LOTAP0059	2
Wepa-Apothekenbedarf	Water	LOTAP0069	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Wepa	Water	LOTAP0072	2
Sophien AM/Wepa	Water	LOTAP0077	1
Wepa	Water	LOTAP0013	1
Wepa	Water	LOTAP0079	2
wepa	Water	LOTAP0024	1
Caelo	Water	LOTAP0039	1
WEPA	Water	LOTAP0078	1
wepa	Water	LOTAP0031	1
Wepa	Water	LOTAP0046	1
Fagron	Water	T14-031432	1
WEPA	Water	LOTAP0051	1

- 88 869 spectra from 501 *Apo-Ident* customers from a total of 26 341 batches from a further 331 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Water and aqueous solutions* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Water and aqueous solutions* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	690	0	34 854
Type B	0	1572	3	91 511
Type C	0	2343	140	88 572

The substance/substance group *Water and aqueous solutions* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.1304 %)
Type B	100.0000 % (> 99.9775 %)	99.8095 % (> 99.6190 %)
Type C	100.0000 % (> 99.6706 %)	94.3617 % (> 94.2408 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested

reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30869	31342	5.49	20.73
31057	32354	2.87	34.55
31240	31463	3.38	41.28
31342	31342	0.00	25.79
31463	31463	0.00	43.52
31624	31718	1.93	38.88
31676	32354	4.01	33.62
31705	31718	9.70	34.18
31718	31718	0.00	38.39
31735	31735	0.00	42.72
31809	31718	4.83	40.81
31829	31463	4.47	40.09
31846	31846	0.00	42.13
31959	31342	5.75	23.64
32308	31718	9.74	35.26
32354	32354	0.00	36.33
32732	31846	2.95	41.82
32831	32831	0.00	36.41

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Wines**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30641-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Wines; Chest elixir; Condurango wine; Elixir e succo liquiritiae; Fortified wine; Gel basis Audor®; Liquiritiae succi elixir; Liquiritiae succus depur. Solut.; Liquorice juice, purified; Succus liquiritiae depur. solut.; Vinum condurangum

Special notes

When selecting the *Wines* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Wines*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Audor	Gel basis Audor®	AP111500	30641	40	20101104
Caelo	Condurango wine	11239009	30901	50	20110909
Audor	Gel basis Audor®	AP103113G	31108	40	1404231

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Fortified wine	13208302	31127	40	1403400
Caelo	Chest elixir	13003109	31213	40	1405332
Caelo	Fortified wine	13208306	31222	70	20130822
Caelo	Fortified wine	13397204	31466	40	20131216
Caelo	Condurango wine	14345603	31581	40	1505325
Caelo	Chest elixir	14182707	31615	40	20140709
Caelo	Liquorice juice,...	14203111	31652	40	20140711
Caelo	Liquorice juice,...	162291	32700	40	20160824*
Audor	Gel basis Audor®	AP53801G	32724	40	20151012
Audor	Gel basis Audor®	AP63101G	32751	40	20160906
Caelo	Liquorice juice,...	171695	33330	40	20170621*
Caelo	Chest elixir	17013603	33619	40	20170125

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 640 spectra of 15 reference samples from the substance/substance group *Wines*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 15 different batches.
- 34 904 spectra from a total of 825 batches from further 207 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 721 spectra from 21 reference samples of the substance/substance group *Wines*.
- These include spectra of independent samples from 14 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Audor Pharma	Gel basis Audor®	zu erfassen	30172	40
Caelo	Fortified wine	10208202	30186	31
Caelo	Liquorice juice, purified	11293204	30415	40
Caelo	Chest elixir	12041501	30504	40
Audor Pharma	Gel basis Audor®	AP103081G	30607	40
Caelo	Condurango wine	11239009	30901 [†]	10
Caelo	Fortified wine	13397204	31466 [†]	20
Caelo	Condurango wine	14345603	31581 [†]	20
Caelo	Chest elixir	14182707	31615 [†]	20
Caelo	Liquorice juice, purified	14203111	31652 [†]	20
Caelo	Liquorice juice, purified	161241	32200	40
Caelo	Fortified wine	161795	32336	40
Audor Pharma	Gel basis Audor®	AP63101G	32752	40
Caelo	Chest elixir	15379113	32956	40
Caelo	Condurango wine	15341903	32957	40
Caelo	Liquorice juice, purified	16316201	33463	40
Audor Pharma	Gel basis Audor®	AP72301G	33626	40
Caelo	Fortified wine	171984	33712	40
Caelo	Condurango wine	15341901	33836	40
Caelo	Condurango wine	180781	34074	40
Caelo	Chest elixir	181495	34228	40

- 93 153 spectra from a total of 2147 batches from further 432 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 123 spectra from 40 *Apo-Ident* customers from 57 batches from the substance/substance group *Wines*.
- These include spectra of independent samples from 54 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Fortified wine	15175501	2
Caelo	Fortified wine	12059809	2
phönix	Fortified wine	16179508	1
Caelo	Fortified wine	16179505	1
Caelo	Fortified wine	12059810	1
Audor Pharma	Gel basis Audor®	Ap103081G	1
Audor Pharma	Gel basis Audor®	AP100852G	4
Audor Pharma	Gel basis Audor®	AP100852	1
Caelo	Chest elixir	3031506	1
Caelo	Chest elixir	12091406	1
Caelo	Chest elixir	18101306	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chest elixir	19071402	1
Caelo	Condurango wine	14345604	1
Audor Pharma	Gel basis Audor®	102464	4
Audor Pharma	Gel basis Audor®	AP102464	1
Audor Pharma	Gel basis Audor®	AP102464G	4
Audor Pharma	Gel basis Audor®		12
Audor Pharma	Gel basis Audor®	AP103993G	3
Audor Pharma	Gel basis Audor®	AP53801G	3
Audor Pharma	Gel basis Audor®	AP72301G	2
Audor Pharma	Gel basis Audor®	ap103113g	1
Fagron	Gel basis Audor®	AP103993G	1
Caelo	Gel basis Audor®	AP53801G	3
Audor Pharma	Gel basis Audor®	AP80601G	2
Caelo	Gel basis Audor®	18041301	1
Anzag	Gel basis Audor®	103081G	1
Großhandel	Gel basis Audor®	AP100852G	1
Audor Pharma	Gel basis Audor®	AP103081	1
Audor Pharma	Gel basis Audor®	AP103081G	2
Audor Pharma	Gel basis Audor®	AP103113G	1
Breko	Fortified wine	201701v12	1
Noweda	Fortified wine	16179509	1
Breko	Fortified wine	201701v13	1
Fiebig	Chest elixir	11019106	1
Caelo	Liquorice juice, purified	13067606	3
Caelo	Liquorice juice, purified	13133801	5
Caelo	Liquorice juice, purified	12116704	5
Caelo	Liquorice juice, purified	13133803	2
Caelo	Liquorice juice, purified	14002602	4
Caelo	Liquorice juice, purified	14203104	2
Caelo	Liquorice juice, purified	14002604	1
Caelo	Liquorice juice, purified	14203107	1
Caelo	Liquorice juice, purified	14203105	3
Caelo	Liquorice juice, purified	14203108	1
Caelo	Liquorice juice, purified	15155302	1
Fagron	Liquorice juice, purified	15155302	2
Caelo	Liquorice juice, purified	16124103	4
Caelo	Liquorice juice, purified	15144203	1
Jenne	Liquorice juice, purified	16124103	1
Caelo	Liquorice juice, purified	16124108	1
Caelo	Liquorice juice, purified	17169508	1
Caelo	Liquorice juice, purified	11293204	3
Caelo	Liquorice juice, purified	16316202	1
Caelo	Liquorice juice, purified	11293208	4
VDL	Liquorice juice, purified	11293210	1
Caelo	Liquorice juice, purified	12061203	1
VDL	Liquorice juice, purified	12116703	2
Caelo	Liquorice juice, purified	12116703	2
Caelo	Liquorice juice, purified	12116709	1
Caelo	Liquorice juice, purified	13067602	2
Caelo	Liquorice juice, purified	13067604	1
Caelo	Liquorice juice, purified	1367604	1
Caelo	Liquorice juice, purified	1367605	1

- 91 229 spectra from 501 *Apo-Ident* customers from a total of 26 952 batches from a further 329 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Wines* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Wines* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	640	0	34 904
Type B	0	711	10	92 786
Type C	0	106	17	91 182

The substance/substance group *Wines* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.0625 %)
Type B	100.0000 % (> 99.9770 %)	98.6130 % (> 98.1969 %)
Type C	100.0000 % (> 99.6588 %)	86.1789 % (> 83.7398 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30641	31108	6.58	24.80
30901	31581	10.46	25.41
31108	31108	0.00	30.66
31127	31127	0.00	23.95
31213	31213	0.00	36.64
31222	31127	2.24	23.48
31466	31127	3.34	22.82
31581	31581	0.00	30.45
31615	31213	8.95	32.07
31652	33330	10.81	34.70
32700	32700	0.00	46.65
32724	31127	9.49	28.14
32751	31127	9.15	22.40
33330	33330	0.00	36.24
33619	31213	6.12	34.84

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Wood oils
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30710-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Wood oils; Cypress oil; Juniper wood oil; Juniperi ligni etheroleum; Oleum cupressi; Oleum Juniperi e Ligno pro uso externo; Oleum terebinthinae rectificat.; Turpentine oil

Special notes

When selecting the *Wood oils* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Wood oils*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Juniper wood oil	11363504	30710	40	20111221
Caelo	Juniper wood oil	11363507	30860	40	20111221

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Juniper wood oil	13196208	31201	40	1405102
Caelo	Turpentine oil	12383210	31214	40	1407275
Euro OTC	Turpentine oil	1312026-01	31531	40	20140113
Euro OTC	Turpentine oil	1511012	31815	40	20151214*
Caelo	Cypress oil	13349109	31837	40	1606484
Caelo	Cypress oil	16130501	32183	40	9464
Euro OTC	Turpentine oil	1703005	33232	40	20170327*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 360 spectra of 9 reference samples from the substance/substance group *Wood oils*. These samples are listed above in the section *calibration samples*. The reference samples originate from 9 different batches.
- 35 184 spectra from a total of 831 batches from further 209 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this *chemometric model*, the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in *Appendix A*. The samples in *Appendix A* were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 475 spectra from 14 reference samples of the substance/substance group *Wood oils*.
- These include spectra of independent samples from 10 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Turpentine oil	10016709	30188	43
Euro OTC	Turpentine oil	1312026-01	31531 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Turpentine oil	1511012	31815 [†]	20
Caelo	Cypress oil	13349109	31837 [†]	20
Caelo	Turpentine oil	161997	32317	40
Caelo	Turpentine oil	171509	33310	40
Caelo	Juniper wood oil	16055410	33506	40
Caelo	Cypress oil	17017101	33602	40
Caelo	Cypress oil	17017101	33603	40
Caelo	Cypress oil	173150	33905	40
Caelo	Juniper wood oil	171964	33914	40
Caelo	Turpentine oil	181108	34060	40
Euro OTC	Turpentine oil	1806017	34188	40
Caelo	Cypress oil	16130505	not required	12

- 93 399 spectra from a total of 2154 batches from further 434 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 86 spectra from 57 *Apo-Ident* customers from 45 batches from the substance/substance group *Wood oils*.
- These include spectra of independent samples from 40 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Juniper wood oil	14008301	1
Caelo	Turpentine oil	15084003	2
Caelo	Turpentine oil	15185101	4
Caelo	Turpentine oil	15185106	3
Caelo	Turpentine oil	11127908	2
Euro OTC	Turpentine oil	711Q-03557	2
Euro OTC	Turpentine oil	1210016-01	4
Caelo	Turpentine oil	12383210	3
Caelo	Turpentine oil	13106701	7
Phönix	Turpentine oil	13378804	1
Caelo	Turpentine oil	13378804	2
Caelo	Turpentine oil	14012601	2
Euro OTC	Turpentine oil	1407037-01	1
Caelo	Turpentine oil	14242703	1
Caelo	Turpentine oil	1424703	1
Euro OTC	Turpentine oil	1503007-01	1
Schiele	Turpentine oil	55379	1
Sanacorp WE:12.01.2017 ...	Turpentine oil	16199701	1
Caelo	Turpentine oil	16199172	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Turpentine oil	161990709	1
Caelo	Turpentine oil	16199701	1
Caelo	Turpentine oil	16199705	2
Caelo	Turpentine oil	16199712	1
Euro OTC	Turpentine oil	981772	2
Caelo	Turpentine oil	412M-02022	1
Caelo	Turpentine oil	10296303	1
Caelo	Turpentine oil	11127905	1
Caelo	Turpentine oil	11358903	4
Caelo	Turpentine oil	11358910	8
Caelo	Turpentine oil	12383201	1
Sanacorp	Turpentine oil	12383210	1
Sanacorp	Turpentine oil	13106701	1
AHC 20.07 13 250ml	Turpentine oil	13106701	1
Caelo	Turpentine oil	14071201	1
Caelo	Turpentine oil	19061310	1
Caelo	Turpentine oil	24544454	1
Caelo	Juniper wood oil	14008306	1
Caelo	Juniper wood oil	14008303	2
Caelo	Juniper wood oil	15149511	1
Caelo	Juniper wood oil	11379007	1
Caelo	Juniper wood oil	12342604	1
Caelo	Juniper wood oil	14008307	3
Caelo	Juniper wood oil	15149502	1
Taoasis GmbH	Juniper wood oil	604898-109389	1
Taoasis GmbH	Juniper wood oil	736525-110967	1
Caelo	Juniper wood oil	11363507	1
Caelo	Juniper wood oil		1
Caelo	Cypress oil		1
Caelo	Cypress oil	16130501	1
Caelo	Cypress oil	17017101	1

- 91 266 spectra from 501 *Apo-Ident* customers from a total of 26 968 batches from a further 331 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Wood oils* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Wood oils* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	360	0	35 184
Type B	0	451	24	93 399
Type C	0	76	10	91 266

The substance/substance group *Wood oils* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.3333 %)
Type B	100.0000 % (> 99.9761 %)	94.9474 % (> 94.3158 %)
Type C	100.0000 % (> 99.6587 %)	88.3721 % (> 84.8837 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30710	31201	3.20	32.23
30860	31201	4.91	31.93
31201	31201	0.00	32.89
31214	31214	0.00	29.50
31531	31214	2.76	30.89
31815	31815	0.00	26.57
31837	31837	0.00	26.56
32183	32183	0.00	26.67
33232	33232	0.00	30.78

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Wool alcohol cream DAB/SR
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30425-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Wool alcohol cream DAB/SR; Lanae alcoholium unguentum aquosum DAB/SR; Unguentum alcoholium lanae aquosum DAB/SR

Special notes

When selecting the *Wool alcohol cream DAB/SR* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Wool alcohol cream DAB/SR*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Wool alcohol cre...	13272603	31229	70	20130829
Caelo	Wool alcohol cre...	13369401	31248	40	20131111
Bombastus	Wool alcohol cre...	284945	31271	40	1406054
Caelo	Wool alcohol cre...	13369401	31287	40	1406055

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Wool alcohol cre...	14031601	31447	40	1512182
Caelo	Wool alcohol cre...	13394101	31455	30	1509376
Caelo	Wool alcohol cre...	14126801	31481	40	20140415
Caelo	Wool alcohol cre...	153608	31923	40	20151029*
Caelo	Wool alcohol cre...	153609	31926	40	20151103*
Caelo	Wool alcohol cre...	153611	31928	40	20151106*
Caelo	Wool alcohol cre...	153613	31930	40	20151106*
Caelo	Wool alcohol cre...	153615	31933	40	20151029*
Caelo	Wool alcohol cre...	153616	31934	40	20151029*
Caelo	Wool alcohol cre...	153617	31937	40	20151103*
Caelo	Wool alcohol cre...	160901	32160	40	20160414*
Caelo	Wool alcohol cre...	170394	33115	40	20170214*
Caelo	Wool alcohol cre...	170396	33140	40	20170214*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 700 spectra of 17 reference samples from the substance/substance group *Wool alcohol cream DAB/SR*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 16 different batches.
- 34 844 spectra from a total of 824 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 4140 spectra from 103 reference samples of the substance/substance group *Wool alcohol cream DAB/SR*.
- These include spectra of independent samples from 79 batches, of which no spectra have been used for the generation of the database.

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Wool alcohol cream DAB/SR	11140801	30317	40
Caelo	Wool alcohol cream DAB/SR	11171701	30405	40
Bombastus	Wool alcohol cream DAB/SR	00263268	30425	40
Caelo	Wool alcohol cream DAB/SR	12066203	30501	40
Caelo	Wool alcohol cream DAB/SR	12066202	30614	40
Bombastus	Wool alcohol cream DAB/SR	00270036	30625	40
Bombastus	Wool alcohol cream DAB/SR	00265302	30627	40
zu erfassen	Wool alcohol cream DAB/SR	zu erfassen	30649	40
Caelo	Wool alcohol cream DAB/SR	12205402	30655	40
Caelo	Wool alcohol cream DAB/SR	12036401	30656	40
apomix	Wool alcohol cream DAB/SR	39A100412	30676	40
Caelo	Wool alcohol cream DAB/SR	153608	31923 [†]	20
Caelo	Wool alcohol cream DAB/SR	153608	31924	60
Caelo	Wool alcohol cream DAB/SR	153608	31925	60
Caelo	Wool alcohol cream DAB/SR	153609	31926 [†]	20
Caelo	Wool alcohol cream DAB/SR	153610	31927	60
Caelo	Wool alcohol cream DAB/SR	153611	31928 [†]	20
Caelo	Wool alcohol cream DAB/SR	153612	31929	60
Caelo	Wool alcohol cream DAB/SR	153613	31930 [†]	20
Caelo	Wool alcohol cream DAB/SR	153614	31931	60
Caelo	Wool alcohol cream DAB/SR	153615	31932	60
Caelo	Wool alcohol cream DAB/SR	153615	31933 [†]	20
Caelo	Wool alcohol cream DAB/SR	153616	31934 [†]	20
Caelo	Wool alcohol cream DAB/SR	153616	31935	60
Caelo	Wool alcohol cream DAB/SR	153617	31936	60
Caelo	Wool alcohol cream DAB/SR	153617	31937 [†]	20
Caelo	Wool alcohol cream DAB/SR	160078	32157	40
Caelo	Wool alcohol cream DAB/SR	160899	32158	40
Caelo	Wool alcohol cream DAB/SR	160900	32159	40
Caelo	Wool alcohol cream DAB/SR	160902	32161	40
Caelo	Wool alcohol cream DAB/SR	160903	32162	40
Caelo	Wool alcohol cream DAB/SR	160904	32163	40
Caelo	Wool alcohol cream DAB/SR	160905	32164	40
Caelo	Wool alcohol cream DAB/SR	160906	32165	40
Caelo	Wool alcohol cream DAB/SR	160907	32166	40
Caelo	Wool alcohol cream DAB/SR	161122	32167	40
Caelo	Wool alcohol cream DAB/SR	160143	32295	40
Caelo	Wool alcohol cream DAB/SR	161208	32296	40
Caelo	Wool alcohol cream DAB/SR	161578	32361	40
Caelo	Wool alcohol cream DAB/SR	161578	32362	40
Caelo	Wool alcohol cream DAB/SR	161577	32363	40
Caelo	Wool alcohol cream DAB/SR	161579	32364	40
Caelo	Wool alcohol cream DAB/SR	161580	32365	40
Caelo	Wool alcohol cream DAB/SR	161580	32366	40
Caelo	Wool alcohol cream DAB/SR	161581	32367	40
Caelo	Wool alcohol cream DAB/SR	161577	32379	40
Caelo	Wool alcohol cream DAB/SR	161776	32380	40
Caelo	Wool alcohol cream DAB/SR	161774	32381	40
Caelo	Wool alcohol cream DAB/SR	161778	32513	40
Caelo	Wool alcohol cream DAB/SR	162490	32514	40
Caelo	Wool alcohol cream DAB/SR	161578	32630	40
Caelo	Wool alcohol cream DAB/SR	161579	32631	40
Caelo	Wool alcohol cream DAB/SR	161992	32632	40
Caelo	Wool alcohol cream DAB/SR	161992	32633	40
Caelo	Wool alcohol cream DAB/SR	161992	32634	40
Caelo	Wool alcohol cream DAB/SR	161993	32635	40
Caelo	Wool alcohol cream DAB/SR	161993	32636	40
Caelo	Wool alcohol cream DAB/SR	161993	32637	40
Caelo	Wool alcohol cream DAB/SR	161994	32638	40
Caelo	Wool alcohol cream DAB/SR	161994	32639	40

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Wool alcohol cream DAB/SR	161994	32640	40
Caelo	Wool alcohol cream DAB/SR	161989	32641	40
Caelo	Wool alcohol cream DAB/SR	163129	32848	40
Caelo	Wool alcohol cream DAB/SR	163131	32849	40
Caelo	Wool alcohol cream DAB/SR	162779	32850	40
Caelo	Wool alcohol cream DAB/SR	162782	32851	40
Caelo	Wool alcohol cream DAB/SR	162783	32852	40
Caelo	Wool alcohol cream DAB/SR	162781	32853	40
Caelo	Wool alcohol cream DAB/SR	162823	32854	40
Caelo	Wool alcohol cream DAB/SR	163130	32887	40
Caelo	Wool alcohol cream DAB/SR	170392	33114	40
Caelo	Wool alcohol cream DAB/SR	170395	33116	40
Caelo	Wool alcohol cream DAB/SR	170398	33117	40
Caelo	Wool alcohol cream DAB/SR	170399	33118	40
Caelo	Wool alcohol cream DAB/SR	170400	33119	40
Caelo	Wool alcohol cream DAB/SR	170401	33120	40
Caelo	Wool alcohol cream DAB/SR	170397	33121	40
Caelo	Wool alcohol cream DAB/SR	163407	33141	40
Caelo	Wool alcohol cream DAB/SR	163409	33142	40
Caelo	Wool alcohol cream DAB/SR	163410	33143	40
Caelo	Wool alcohol cream DAB/SR	163441	33186	40
Caelo	Wool alcohol cream DAB/SR	170941	33372	40
Caelo	Wool alcohol cream DAB/SR	171028	33405	40
Caelo	Wool alcohol cream DAB/SR	171027	33406	40
Caelo	Wool alcohol cream DAB/SR	171026	33407	40
Caelo	Wool alcohol cream DAB/SR	171487	33518	40
Caelo	Wool alcohol cream DAB/SR	171490	33519	40
Caelo	Wool alcohol cream DAB/SR	171492	33520	40
Caelo	Wool alcohol cream DAB/SR	171876	33521	40
Caelo	Wool alcohol cream DAB/SR	171877	33522	40
Caelo	Wool alcohol cream DAB/SR	171878	33523	40
Caelo	Wool alcohol cream DAB/SR	171879	33524	40
Caelo	Wool alcohol cream DAB/SR	171880	33525	40
Caelo	Wool alcohol cream DAB/SR	171783	33549	40
Caelo	Wool alcohol cream DAB/SR	171782	33758	40
Caelo	Wool alcohol cream DAB/SR	172844	33762	40
Caelo	Wool alcohol cream DAB/SR	172845	33763	40
Caelo	Wool alcohol cream DAB/SR	172846	33764	40
Caelo	Wool alcohol cream DAB/SR	172847	33765	40
Caelo	Wool alcohol cream DAB/SR	172848	33766	40
Caelo	Wool alcohol cream DAB/SR	172287	33770	40
Caelo	Wool alcohol cream DAB/SR	180234	34156	40
Caelo	Wool alcohol cream DAB/SR	180863	34157	40

- 89 734 spectra from a total of 2081 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 3217 spectra from 391 *Apo-Ident* customers from 1067 batches from the substance/substance group *Wool alcohol cream DAB/SR*.
- These include spectra of independent samples from 1057 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	116282301	1
Caelo	Wool alcohol cream DAB/SR	15360801	4
Caelo	Wool alcohol cream DAB/SR	13145201	7
Caelo	Wool alcohol cream DAB/SR	15360802	5
Caelo	Wool alcohol cream DAB/SR	15361101	4
Caelo	Wool alcohol cream DAB/SR	15361201	10
Caelo	Wool alcohol cream DAB/SR	15361402	11
Caelo	Wool alcohol cream DAB/SR	15361501	8
Caelo	Wool alcohol cream DAB/SR	15361601	3
Caelo	Wool alcohol cream DAB/SR	16007701	10
Caelo	Wool alcohol cream DAB/SR	16007901	14
Phoenix	Wool alcohol cream DAB/SR	16007701	1
Caelo	Wool alcohol cream DAB/SR	16008103	21
Caelo	Wool alcohol cream DAB/SR	16008201	5
Caelo	Wool alcohol cream DAB/SR	16008301	6
Caelo	Wool alcohol cream DAB/SR	16008401	1
Caesar & Loretz/Noweda	Wool alcohol cream DAB/SR	16008401	1
Sonnen-Apotheke Wefensl...	Wool alcohol cream DAB/SR	16090101	1
Caelo	Wool alcohol cream DAB/SR	16089902	6
WE:Sanacorp26.09.2016 ...	Wool alcohol cream DAB/SR	16090002	1
Caelo	Wool alcohol cream DAB/SR	16090101	2
Caelo	Wool alcohol cream DAB/SR	16090201	1
Caelo	Wool alcohol cream DAB/SR	16090301	3
Caelo	Wool alcohol cream DAB/SR	16090401	8
Caelo	Wool alcohol cream DAB/SR	16090501	21
Caelo	Wool alcohol cream DAB/SR	16090702	15
29.06.2016 Phoenix 3,77EUR	Wool alcohol cream DAB/SR	16090702	1
Caelo	Wool alcohol cream DAB/SR	160912ECA	1
Caelo	Wool alcohol cream DAB/SR	16112201	8
Caelo	Wool alcohol cream DAB/SR	16157801	9
Sanacorp/AHD	Wool alcohol cream DAB/SR	16157801	1
Caelo	Wool alcohol cream DAB/SR	16157802	6
Caelo	Wool alcohol cream DAB/SR	16157901	14
Caesar & Loretz GmbH Hi...	Wool alcohol cream DAB/SR	16157901	1
Bombastus	Wool alcohol cream DAB/SR	16158001	2
Caelo	Wool alcohol cream DAB/SR	16158001	10
Caelo	Wool alcohol cream DAB/SR	16158101	7
Caelo	Wool alcohol cream DAB/SR	16158102	3
Caelo	Wool alcohol cream DAB/SR	16199003	1
Caelo	Wool alcohol cream DAB/SR	16199004	9
Celo	Wool alcohol cream DAB/SR	16199004	1
Euro OTC	Wool alcohol cream DAB/SR	16199004	1
Caelo	Wool alcohol cream DAB/SR	16199101	3
Caelo	Wool alcohol cream DAB/SR	16199201	5
Bombastus	Wool alcohol cream DAB/SR	16199201	1
SanacorpWE:11.01.2017 ...	Wool alcohol cream DAB/SR	16199201	1
Caelo	Wool alcohol cream DAB/SR	16199301	2
Caesar & Loretz/Noweda	Wool alcohol cream DAB/SR	16199401	1
Caelo	Wool alcohol cream DAB/SR	16277902	6
Caelo	Wool alcohol cream DAB/SR	16277903	5
Caelo	Wool alcohol cream DAB/SR	16278201	2
Caelo	Wool alcohol cream DAB/SR	16278301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	16282301	19
Caelo	Wool alcohol cream DAB/SR	163123101	1
Caelo	Wool alcohol cream DAB/SR	16312901	4
Caelo	Wool alcohol cream DAB/SR	16312902	6
Caelo	Wool alcohol cream DAB/SR	16313002	3
Caelo	Wool alcohol cream DAB/SR	16313101	11
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	16313101	1
Caelo	Wool alcohol cream DAB/SR	16340703	2
Caelo	Wool alcohol cream DAB/SR	16340701	2
Noweda	Wool alcohol cream DAB/SR	16313101	1
Caelo	Wool alcohol cream DAB/SR	16340801	17
Noweda	Wool alcohol cream DAB/SR	16340801	1
Caelo	Wool alcohol cream DAB/SR	16340902	2
Caelo	Wool alcohol cream DAB/SR	16341001	3
Caelo	Wool alcohol cream DAB/SR	1703401	1
Caelo	Wool alcohol cream DAB/SR	170390301	1
Caelo	Wool alcohol cream DAB/SR	17039201	4
Caelo	Wool alcohol cream DAB/SR	17039301	10
Finteler Apotheke	Wool alcohol cream DAB/SR	17039301	1
Caelo	Wool alcohol cream DAB/SR	17039401	9
Wepa/ Phoenix	Wool alcohol cream DAB/SR	17039401	1
Caelo	Wool alcohol cream DAB/SR	17039701	3
Caelo	Wool alcohol cream DAB/SR	17039801	2
Caelo	Wool alcohol cream DAB/SR	17039901	1
Caelo	Wool alcohol cream DAB/SR	17040001	1
Caelo	Wool alcohol cream DAB/SR	17094001	2
Fagron	Wool alcohol cream DAB/SR	17040101	1
Caelo	Wool alcohol cream DAB/SR	17102602	1
Sanacorp WE: 04.09.17 ...	Wool alcohol cream DAB/SR	17102603	1
Noweda	Wool alcohol cream DAB/SR	17102701	1
Caelo	Wool alcohol cream DAB/SR	17102701	10
Caelo	Wool alcohol cream DAB/SR	17102603	2
Caesar&Loretz	Wool alcohol cream DAB/SR	17040101	1
Caelo	Wool alcohol cream DAB/SR	17128502	2
Caelo	Wool alcohol cream DAB/SR	17148901	2
Caelo	Wool alcohol cream DAB/SR	17149002	3
Caelo	Wool alcohol cream DAB/SR	17155901	5
Caelo	Wool alcohol cream DAB/SR	17149101	1
Caelo	Wool alcohol cream DAB/SR	17148701	2
Noweda	Wool alcohol cream DAB/SR	17156001	1
Caelo	Wool alcohol cream DAB/SR	17156001	5
Caelo	Wool alcohol cream DAB/SR	17187701	5
Caelo	Wool alcohol cream DAB/SR	17187801	4
Noweda	Wool alcohol cream DAB/SR	17187801	1
Caelo	Wool alcohol cream DAB/SR	17187902	1
Caelo	Wool alcohol cream DAB/SR	17284301	2
Sanacorp W.E: 24.10.18 ...	Wool alcohol cream DAB/SR	17188001	1
Caelo	Wool alcohol cream DAB/SR	17284502	1
Caelo	Wool alcohol cream DAB/SR	17284701	1
Caelo	Wool alcohol cream DAB/SR	18012201	3
Caelo	Wool alcohol cream DAB/SR	17284802	2
Caelo	Wool alcohol cream DAB/SR	18096001	1
Caelo	Wool alcohol cream DAB/SR	20170419ECA/16312902	1
Caelo	Wool alcohol cream DAB/SR	19199004	1
Caelo	Wool alcohol cream DAB/SR	2716Q-01371	1
Caelo	Wool alcohol cream DAB/SR	18012601	1
Beiersdorf/Alliance	Wool alcohol cream DAB/SR	51630108WA	1
Caelo	Wool alcohol cream DAB/SR	18012701	1
Caelo	Wool alcohol cream DAB/SR	14362201	18
Caelo	Wool alcohol cream DAB/SR	14259701	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	14259601	11
Caelo	Wool alcohol cream DAB/SR	1015E-01371	1
Caelo	Wool alcohol cream DAB/SR	14362501	8
Caelo	Wool alcohol cream DAB/SR	14295601	6
Caelo	Wool alcohol cream DAB/SR	15017502	5
Caelo	Wool alcohol cream DAB/SR	15017601	21
Caelo	Wool alcohol cream DAB/SR	15017901	7
Caelo	Wool alcohol cream DAB/SR	15017801	3
Caelo	Wool alcohol cream DAB/SR	15073601	17
Caelo	Wool alcohol cream DAB/SR	15073701	3
Caelo	Wool alcohol cream DAB/SR	15073801	2
Caelo	Wool alcohol cream DAB/SR	15149801	12
Sanacorp /Caelo 5.12.15	Wool alcohol cream DAB/SR	15149801	2
Caelo/ Noweda	Wool alcohol cream DAB/SR	15150001	2
Caelo	Wool alcohol cream DAB/SR	15150001	6
Caelo	Wool alcohol cream DAB/SR	15150002	2
Caelo	Wool alcohol cream DAB/SR	15150101	7
Caelo	Wool alcohol cream DAB/SR	15226301	8
Gehe	Wool alcohol cream DAB/SR	15226301	1
Caelo	Wool alcohol cream DAB/SR	15259701	19
Sanacorp /Caelo 7.12.15	Wool alcohol cream DAB/SR	15259701	1
Caelo	Wool alcohol cream DAB/SR	15259801	6
Caelo	Wool alcohol cream DAB/SR	15259802	1
Caelo	Wool alcohol cream DAB/SR	15259901	5
Caelo	Wool alcohol cream DAB/SR	15260001	16
Caelo/Noweda	Wool alcohol cream DAB/SR	15259901	1
Caelo	Wool alcohol cream DAB/SR	15260102	20
Caelo	Wool alcohol cream DAB/SR	15266001	6
Caelo	Wool alcohol cream DAB/SR	15266101	2
Bombastus	Wool alcohol cream DAB/SR	15266001	1
Caelo	Wool alcohol cream DAB/SR	15266201	7
Caelo	Wool alcohol cream DAB/SR	1526601	1
Sanacorp	Wool alcohol cream DAB/SR	15266201	2
Caesar&Loretz/Caesar&Lo...	Wool alcohol cream DAB/SR	15266301	1
Caelo	Wool alcohol cream DAB/SR	15266301	5
Caelo	Wool alcohol cream DAB/SR	15266401	4
Caelo	Wool alcohol cream DAB/SR	15266501	6
unisan	Wool alcohol cream DAB/SR	15266501	1
Caelo	Wool alcohol cream DAB/SR	1529701	1
Caelo	Wool alcohol cream DAB/SR	15360901	15
Caelo	Wool alcohol cream DAB/SR	15360902	4
Caelo	Wool alcohol cream DAB/SR	15361001	1
Caesar & Loretz GmbH Hi...	Wool alcohol cream DAB/SR	15361001	2
Caelo	Wool alcohol cream DAB/SR	15361301	6
WE:05.04.2016 Sanacorp	Wool alcohol cream DAB/SR	15361501	1
Caelo	Wool alcohol cream DAB/SR	153671601	1
WE:17.02.2016	Wool alcohol cream DAB/SR	15361601	1
Caelo	Wool alcohol cream DAB/SR	16007702	3
Caelo	Wool alcohol cream DAB/SR	16007801	1
Caelo	Wool alcohol cream DAB/SR	16007802	2
Noweda/ Caesar & Loretz...	Wool alcohol cream DAB/SR	16007901	1
Fagron	Wool alcohol cream DAB/SR	16007901	1
Fiebig	Wool alcohol cream DAB/SR	16008201	1
Caelo	Wool alcohol cream DAB/SR	1608001	1
Sanacorp	Wool alcohol cream DAB/SR	16090703	1
Sanacorp	Wool alcohol cream DAB/SR	27/AKD4	1
Caelo	Wool alcohol cream DAB/SR		6
Bombastus	Wool alcohol cream DAB/SR	295460	2
Caelo	Wool alcohol cream DAB/SR	1118956	1
Caelo	Wool alcohol cream DAB/SR	269213261801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	614E1371	1
Caelo	Wool alcohol cream DAB/SR	80616-10	1
Caelo	Wool alcohol cream DAB/SR	9071301	1
Caelo	Wool alcohol cream DAB/SR	1098277	1
Caelo	Wool alcohol cream DAB/SR	11171501	1
Caelo	Wool alcohol cream DAB/SR	11228001	2
Ichthyol-Gesellschaft	Wool alcohol cream DAB/SR	11288302	1
Caelo	Wool alcohol cream DAB/SR	11356001	2
Caelo	Wool alcohol cream DAB/SR	11356101	3
Caelo	Wool alcohol cream DAB/SR	11356301	1
Caelo	Wool alcohol cream DAB/SR	11356201	1
Caelo	Wool alcohol cream DAB/SR	11356601	1
Caelo	Wool alcohol cream DAB/SR	12000201	2
Caelo	Wool alcohol cream DAB/SR	12001101	3
Caelo	Wool alcohol cream DAB/SR	12000301	1
Caelo	Wool alcohol cream DAB/SR	12001301	1
Caelo	Wool alcohol cream DAB/SR	12030701	1
Caelo	Wool alcohol cream DAB/SR	12030601	2
Caelo	Wool alcohol cream DAB/SR	12001201	1
Caelo	Wool alcohol cream DAB/SR	11356501	1
Caelo	Wool alcohol cream DAB/SR	12036601	5
Caelo	Wool alcohol cream DAB/SR	120307010	1
Caelo	Wool alcohol cream DAB/SR	12036501	5
Caelo	Wool alcohol cream DAB/SR	12056202	1
Caelo	Wool alcohol cream DAB/SR	12056301	1
Caelo	Wool alcohol cream DAB/SR	12056401	3
ANZAG am 18.05.2012	Wool alcohol cream DAB/SR	12056402	1
Caelo	Wool alcohol cream DAB/SR	12101304	1
Caelo	Wool alcohol cream DAB/SR	12106001	4
Caelo	Wool alcohol cream DAB/SR	12106102	3
Caelo	Wool alcohol cream DAB/SR	12106601	1
Caelo	Wool alcohol cream DAB/SR	12106401	3
Devesa	Wool alcohol cream DAB/SR	12106601	1
Noweda	Wool alcohol cream DAB/SR	12121901	3
Gehe	Wool alcohol cream DAB/SR	12122001	1
Caelo	Wool alcohol cream DAB/SR	12122001	5
Caelo	Wool alcohol cream DAB/SR	12106701	2
Caelo	Wool alcohol cream DAB/SR	12122002	4
Caelo	Wool alcohol cream DAB/SR	1211901	1
Caelo	Wool alcohol cream DAB/SR	1212202	1
Caelo	Wool alcohol cream DAB/SR	12122101	6
Caelo	Wool alcohol cream DAB/SR	12122201	1
Caelo	Wool alcohol cream DAB/SR	12123501	6
Caelo	Wool alcohol cream DAB/SR	12205401	1
Caelo	Wool alcohol cream DAB/SR	12205402	1
Caelo	Wool alcohol cream DAB/SR	12205501	1
Caelo	Wool alcohol cream DAB/SR	12205701	6
Anzag 04.12.12	Wool alcohol cream DAB/SR	12205402	1
Caelo	Wool alcohol cream DAB/SR	12208902	8
Caelo	Wool alcohol cream DAB/SR	12205801	3
Noweda	Wool alcohol cream DAB/SR	12208902	1
Post-Apotheke	Wool alcohol cream DAB/SR	12209001	2
Caelo	Wool alcohol cream DAB/SR	12209001	2
anzag, 05.11.2012	Wool alcohol cream DAB/SR	12209101	1
Caelo	Wool alcohol cream DAB/SR	12209101	3
Caelo	Wool alcohol cream DAB/SR	12259701	3
Caelo	Wool alcohol cream DAB/SR	12259702	4
Caelo	Wool alcohol cream DAB/SR	12259802	2
Caelo	Wool alcohol cream DAB/SR	12259901	3
Caelo	Wool alcohol cream DAB/SR	12260002	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Wool alcohol cream DAB/SR	12260002	1
Sanacorp	Wool alcohol cream DAB/SR	12260101	1
Caelo	Wool alcohol cream DAB/SR	12260402	15
Caelo	Wool alcohol cream DAB/SR	12260101	2
Anzag 27.12.12 ek:9,46	Wool alcohol cream DAB/SR	12260101	1
Caelo	Wool alcohol cream DAB/SR	12269201	4
Caelo	Wool alcohol cream DAB/SR	12307102	1
Caelo	Wool alcohol cream DAB/SR	12307103	4
phönix 22.01.2013	Wool alcohol cream DAB/SR	12307103	1
Caelo	Wool alcohol cream DAB/SR	12307201	8
Sanacorp	Wool alcohol cream DAB/SR	12307201	1
Caelo	Wool alcohol cream DAB/SR	12307301	1
Caelo	Wool alcohol cream DAB/SR	12307401	7
Gehe	Wool alcohol cream DAB/SR	12307401	1
Caelo	Wool alcohol cream DAB/SR	12307501	4
Sanacorp	Wool alcohol cream DAB/SR	12307501	1
Caelo	Wool alcohol cream DAB/SR	12307601	5
Großhandel	Wool alcohol cream DAB/SR	12307601	1
Caelo	Wool alcohol cream DAB/SR	12307701	6
Caelo	Wool alcohol cream DAB/SR	12307801	3
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	12307901	1
Caelo	Wool alcohol cream DAB/SR	12307901	2
Phoenix 27.02.13	Wool alcohol cream DAB/SR	12308001	1
Caelo	Wool alcohol cream DAB/SR	12308001	7
Caelo	Wool alcohol cream DAB/SR	12351901	5
Fagron	Wool alcohol cream DAB/SR	12370202	1
Caelo	Wool alcohol cream DAB/SR	12370202	2
Caelo	Wool alcohol cream DAB/SR	12370203	1
Caelo	Wool alcohol cream DAB/SR	12370401	12
Caelo	Wool alcohol cream DAB/SR	12370601	9
Phoenix 11.03.13	Wool alcohol cream DAB/SR	1238001	1
Großhandel	Wool alcohol cream DAB/SR	12370601	1
Caelo	Wool alcohol cream DAB/SR	13029502	3
Caelo	Wool alcohol cream DAB/SR	13029503	6
Finze	Wool alcohol cream DAB/SR	13029503	1
Caelo	Wool alcohol cream DAB/SR	13061401	3
AHD 09.04.2013	Wool alcohol cream DAB/SR	13029504	1
Caelo	Wool alcohol cream DAB/SR	13061402	3
Caelo	Wool alcohol cream DAB/SR	13061501	14
Gehe	Wool alcohol cream DAB/SR	13061501	1
Fagron	Wool alcohol cream DAB/SR	13061501	1
Caelo	Wool alcohol cream DAB/SR	13061601	4
Caelo	Wool alcohol cream DAB/SR	13061701	2
Ichthyol/Phönix	Wool alcohol cream DAB/SR	13061801	2
Fagron	Wool alcohol cream DAB/SR	13061701	1
Großhandel	Wool alcohol cream DAB/SR	13061901	1
Caelo	Wool alcohol cream DAB/SR	13061901	7
Caelo	Wool alcohol cream DAB/SR	13062101	13
Sanacorp	Wool alcohol cream DAB/SR	13062101	2
Post-Apotheke	Wool alcohol cream DAB/SR	13062001	1
Caelo	Wool alcohol cream DAB/SR	13062201	3
Gehe 02.07.2013	Wool alcohol cream DAB/SR	13062201	3
Gehe	Wool alcohol cream DAB/SR	13062101	1
Cealo/ Noweda	Wool alcohol cream DAB/SR	13068504	2
Caelo	Wool alcohol cream DAB/SR	13066501	1
Caelo	Wool alcohol cream DAB/SR	13062301	5
Caelo	Wool alcohol cream DAB/SR	13068701	9
Caesar & Loretz/ Noweda	Wool alcohol cream DAB/SR	13145202	1
Caelo	Wool alcohol cream DAB/SR	13145202	6
Caelo	Wool alcohol cream DAB/SR	13145301	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	13145401	4
A.H. 25.06.13 ek: 9,75EUR	Wool alcohol cream DAB/SR	13145501	1
Großhandel	Wool alcohol cream DAB/SR	13145501	1
Caelo	Wool alcohol cream DAB/SR	13145501	5
Caelo	Wool alcohol cream DAB/SR	13145601	14
Anzag	Wool alcohol cream DAB/SR	13145801	2
Caelo	Wool alcohol cream DAB/SR	13145701	3
Caelo	Wool alcohol cream DAB/SR	13145801	7
DEFEKTUR	Wool alcohol cream DAB/SR	13145801	1
Caelo	Wool alcohol cream DAB/SR	13145901	5
Caelo	Wool alcohol cream DAB/SR	13146001	2
Bombastus	Wool alcohol cream DAB/SR	13145901	1
Gehe 26.09.13	Wool alcohol cream DAB/SR	13146001	1
Caesar & Loretz GmbH, 4...	Wool alcohol cream DAB/SR	13171901	1
Caelo	Wool alcohol cream DAB/SR	13190201	1
Großhandel	Wool alcohol cream DAB/SR	13190203	1
Caelo	Wool alcohol cream DAB/SR	13261701	13
Caelo	Wool alcohol cream DAB/SR	13172103	11
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	13261701	1
Caelo	Wool alcohol cream DAB/SR	13261801	13
Caelo	Wool alcohol cream DAB/SR	13279001	6
Caelo	Wool alcohol cream DAB/SR	13279002	2
Gehe	Wool alcohol cream DAB/SR	13279002	1
Gehe 24.10.2013	Wool alcohol cream DAB/SR	13279002	1
Caelo	Wool alcohol cream DAB/SR	13279101	4
Caelo	Wool alcohol cream DAB/SR	13279201	4
Caelo	Wool alcohol cream DAB/SR	13279301	4
Phönix	Wool alcohol cream DAB/SR	13279301	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	13279301	1
Caelo	Wool alcohol cream DAB/SR	13279401	3
Caelo	Wool alcohol cream DAB/SR	13279502	4
Caelo	Wool alcohol cream DAB/SR	13279701	2
Caelo	Wool alcohol cream DAB/SR	13289301	10
Fagron	Wool alcohol cream DAB/SR	13289301	1
Gehe am 03.02.2014	Wool alcohol cream DAB/SR	13289301	1
Caelo	Wool alcohol cream DAB/SR	13289403	7
Gehe am 14.01.2014	Wool alcohol cream DAB/SR	13289403	1
Caelo	Wool alcohol cream DAB/SR	13322801	22
Noweda	Wool alcohol cream DAB/SR	13322801	1
Caelo	Wool alcohol cream DAB/SR	13369401	8
Caelo	Wool alcohol cream DAB/SR	1336901	1
Caelo	Wool alcohol cream DAB/SR	13369601	15
Sanacorp	Wool alcohol cream DAB/SR	13369601	1
Caelo	Wool alcohol cream DAB/SR	13369701	4
Alliance Health 07.01.2...	Wool alcohol cream DAB/SR	13369701	1
Caelo	Wool alcohol cream DAB/SR	13369801	5
Caelo	Wool alcohol cream DAB/SR	13373501	3
Caelo	Wool alcohol cream DAB/SR	13373601	5
Caelo	Wool alcohol cream DAB/SR	13373701	10
Caelo	Wool alcohol cream DAB/SR	13373903	3
Caelo	Wool alcohol cream DAB/SR	13394101	10
Caelo	Wool alcohol cream DAB/SR	13394102	22
Noweda	Wool alcohol cream DAB/SR	14002801	2
Sanacorp	Wool alcohol cream DAB/SR	14002801	1
Caelo	Wool alcohol cream DAB/SR	14002801	11
Caelo	Wool alcohol cream DAB/SR	14045602	1
Caelo	Wool alcohol cream DAB/SR	14033202	7
Caelo	Wool alcohol cream DAB/SR	14045702	3
Caelo	Wool alcohol cream DAB/SR	14046201	5
Caelo	Wool alcohol cream DAB/SR	14126801	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	14046401	5
Caelo	Wool alcohol cream DAB/SR	14164201	18
Caelo	Wool alcohol cream DAB/SR	14162503	4
Phoenix	Wool alcohol cream DAB/SR	14164201	1
Caelo	Wool alcohol cream DAB/SR	14162502	11
Caelo	Wool alcohol cream DAB/SR	14167101	5
Caelo	Wool alcohol cream DAB/SR	14166901	4
Caelo	Wool alcohol cream DAB/SR	14169101	11
Caelo	Wool alcohol cream DAB/SR	14259801	11
Caelo	Wool alcohol cream DAB/SR	14170801	3
Caelo	Wool alcohol cream DAB/SR	1426401	1
Caelo	Wool alcohol cream DAB/SR	1429501	2
Caelo	Wool alcohol cream DAB/SR	14295801	6
Ichthyol Gesellschaft C...	Wool alcohol cream DAB/SR	14362501	1
Caelo	Wool alcohol cream DAB/SR	15149802	2
Phönix	Wool alcohol cream DAB/SR	15259701	2
Celo	Wool alcohol cream DAB/SR	16008301	1
Caelo	Wool alcohol cream DAB/SR	1608992	1
Caelo	Wool alcohol cream DAB/SR	16090601	1
Caelo	Wool alcohol cream DAB/SR	1913Q-01371	1
Caelo	Wool alcohol cream DAB/SR	2x13279201	1
DEFEKTUR	Wool alcohol cream DAB/SR	2013/08/14	1
Caelo	Wool alcohol cream DAB/SR	22/AKD52	1
Post-Apotheke	Wool alcohol cream DAB/SR	5113F-07109	1
Caelo	Wool alcohol cream DAB/SR	Ch12307102	1
Caelo	Wool alcohol cream DAB/SR	E-090413-3	1
Caelo	Wool alcohol cream DAB/SR	PE113	1
Apomix	Wool alcohol cream DAB/SR	3A160113	1
Caelo	Wool alcohol cream DAB/SR	13373801	8
VDL 28.04;Preis:4,34	Wool alcohol cream DAB/SR		1
Caelo	Wool alcohol cream DAB/SR	14045801	13
Caelo	Wool alcohol cream DAB/SR	14246101	5
Caelo	Wool alcohol cream DAB/SR	14246201	4
Caelo	Wool alcohol cream DAB/SR	14295401	5
Bombastus	Wool alcohol cream DAB/SR	284451	3
Bombastus	Wool alcohol cream DAB/SR	284582	1
Bombastus	Wool alcohol cream DAB/SR	285868	2
Bombastus	Wool alcohol cream DAB/SR	289393	1
Bombastus	Wool alcohol cream DAB/SR	289384	1
Bombastus	Wool alcohol cream DAB/SR	289489	2
Bombastus	Wool alcohol cream DAB/SR	288442	2
Bombastus	Wool alcohol cream DAB/SR	287087	2
Bombastus	Wool alcohol cream DAB/SR	291264	1
Caelo	Wool alcohol cream DAB/SR	13029504	1
Caelo	Wool alcohol cream DAB/SR	1314570	1
Sanacorp	Wool alcohol cream DAB/SR	13145601	1
Caelo	Wool alcohol cream DAB/SR	13279201-3433311180	1
Caelo	Wool alcohol cream DAB/SR	13369402	2
Bombastus	Wool alcohol cream DAB/SR	274634	1
Bombastus	Wool alcohol cream DAB/SR	275990	1
Bombastus	Wool alcohol cream DAB/SR	277206	1
Bombastus	Wool alcohol cream DAB/SR	278745	1
Bombastus	Wool alcohol cream DAB/SR	275911	1
Bombastus	Wool alcohol cream DAB/SR	278827	1
Bombastus	Wool alcohol cream DAB/SR	279875	2
Bombastus	Wool alcohol cream DAB/SR	278742	1
Bombastus	Wool alcohol cream DAB/SR	280386	1
Bombastus	Wool alcohol cream DAB/SR	280390	2
Caelo	Wool alcohol cream DAB/SR	14092502	7
Vdl 10.07.14	Wool alcohol cream DAB/SR	14126801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	14167001	4
Caelo	Wool alcohol cream DAB/SR	15017501	3
Sanacorp	Wool alcohol cream DAB/SR	15017601	1
VDL;4,87EUR;12.06.15	Wool alcohol cream DAB/SR	15073601	1
Bombastus	Wool alcohol cream DAB/SR	294270	4
Bombastus	Wool alcohol cream DAB/SR	294330	2
Bombastus	Wool alcohol cream DAB/SR	60080302	2
Bombastus	Wool alcohol cream DAB/SR	267046	1
Bombastus	Wool alcohol cream DAB/SR	269832	1
Anzag	Wool alcohol cream DAB/SR	12122201	1
Caelo	Wool alcohol cream DAB/SR	1206202	1
Anzag	Wool alcohol cream DAB/SR	12269201	1
Anzag	Wool alcohol cream DAB/SR	13029503	1
Caelo	Wool alcohol cream DAB/SR	13171901	7
Caelo	Wool alcohol cream DAB/SR	13068601	8
Gehe	Wool alcohol cream DAB/SR	272253	1
Bombastus	Wool alcohol cream DAB/SR	273631	2
Bombastus	Wool alcohol cream DAB/SR	272259	2
Bombastus	Wool alcohol cream DAB/SR	272263	1
Gehe	Wool alcohol cream DAB/SR	277207	1
Bombastus	Wool alcohol cream DAB/SR	277208	1
Gehe	Wool alcohol cream DAB/SR	282245	1
Bombastus	Wool alcohol cream DAB/SR	285875	1
Bombastus	Wool alcohol cream DAB/SR	284447	1
Gehe	Wool alcohol cream DAB/SR	288442	2
Bombastus	Wool alcohol cream DAB/SR	282685	1
Apomix	Wool alcohol cream DAB/SR	1A020114	1
Caelo	Wool alcohol cream DAB/SR	1071503	1
Pharmachemie	Wool alcohol cream DAB/SR	31215	1
apomix/Noweda	Wool alcohol cream DAB/SR	2A120117	1
pharmachem gmbh co kg,p...	Wool alcohol cream DAB/SR	31215	1
Allianc. Health. / Apomix	Wool alcohol cream DAB/SR	2A120117	1
apomix	Wool alcohol cream DAB/SR	3A270116	1
Caelo	Wool alcohol cream DAB/SR	415A-01371	1
apomix / Anzag	Wool alcohol cream DAB/SR	4A220115	1
Allianc. Health. / Apomix	Wool alcohol cream DAB/SR	5A220217	1
apomix/GEHE23.04.14	Wool alcohol cream DAB/SR	5A200114	1
Apomix/Phönix	Wool alcohol cream DAB/SR	5A220115	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	6A220115	1
apomix	Wool alcohol cream DAB/SR	5A220217	1
Apomix	Wool alcohol cream DAB/SR	6A230217	1
Apomix	Wool alcohol cream DAB/SR	5A280116	1
Apomix /Noweda	Wool alcohol cream DAB/SR	6a230217	1
Gehe/Apomix	Wool alcohol cream DAB/SR	8A290114	1
Caelo	Wool alcohol cream DAB/SR	9041402	1
Apomix Anzag	Wool alcohol cream DAB/SR	9A050215	1
Apomix	Wool alcohol cream DAB/SR	100A021214	1
apomix / Alliance	Wool alcohol cream DAB/SR	100A191213	1
Apomix / Phönix	Wool alcohol cream DAB/SR	101A191213	1
apomix/AHD	Wool alcohol cream DAB/SR	102A091215	1
Apomix/Gehe	Wool alcohol cream DAB/SR	104A091215	1
apomix/Gehe	Wool alcohol cream DAB/SR	104A111214	1
Apomix/Noweda	Wool alcohol cream DAB/SR	105A301214	1
Phönix	Wool alcohol cream DAB/SR	1098283	1
Pharmachem/Phönix	Wool alcohol cream DAB/SR	110614	1
Apomix	Wool alcohol cream DAB/SR	108A301214	1
Pharmachem.Pößneck/Phoe...	Wool alcohol cream DAB/SR	110913	1
Pharmachem GmbH	Wool alcohol cream DAB/SR	110913	1
apomix/Holdermann	Wool alcohol cream DAB/SR	11A270417	1
Caelo	Wool alcohol cream DAB/SR	114163901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix/Noweda	Wool alcohol cream DAB/SR	12A040214	1
Apomix/Noweda	Wool alcohol cream DAB/SR	12A050215	1
Caelo	Wool alcohol cream DAB/SR	13172102	1
Caelo	Wool alcohol cream DAB/SR	13013201	1
Sanacorp	Wool alcohol cream DAB/SR	13172102	1
Caelo	Wool alcohol cream DAB/SR	13272601	1
Caelo	Wool alcohol cream DAB/SR	13280701	4
Caelo	Wool alcohol cream DAB/SR	13281003	4
Caelo	Wool alcohol cream DAB/SR	13369501	2
Großhandel	Wool alcohol cream DAB/SR	13373501	1
Caesar & Loretz GmbH Hi...	Wool alcohol cream DAB/SR	13373701	1
GEHE 24.02.2014	Wool alcohol cream DAB/SR	13373801	1
Caelo	Wool alcohol cream DAB/SR	1339401	1
Caelo	Wool alcohol cream DAB/SR	13394103	2
Caelo	Wool alcohol cream DAB/SR	13394201	7
Gehe	Wool alcohol cream DAB/SR	13394201	1
Caelo	Wool alcohol cream DAB/SR	1373801	1
Apomix Anzag	Wool alcohol cream DAB/SR	13A010617	3
Apomix	Wool alcohol cream DAB/SR	13A050215	1
Apomix/ Gehe	Wool alcohol cream DAB/SR	13A150216	1
apomix	Wool alcohol cream DAB/SR	13A150216	1
Apomix /Alliance	Wool alcohol cream DAB/SR	13A180214	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	14045702	1
Ichthyol Gesellschaft	Wool alcohol cream DAB/SR	14045801	1
Caelo	Wool alcohol cream DAB/SR	14045901	4
Caelo	Wool alcohol cream DAB/SR	14046001	1
Caelo	Wool alcohol cream DAB/SR	14046002	1
Caelo	Wool alcohol cream DAB/SR	14046101	2
Phönix	Wool alcohol cream DAB/SR	14046101	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	14046201	1
Sanacorp	Wool alcohol cream DAB/SR	14046401	1
Caelo	Wool alcohol cream DAB/SR	14046501	4
Caelo	Wool alcohol cream DAB/SR	14046502	3
Caelo	Wool alcohol cream DAB/SR	14092401	3
Caelo	Wool alcohol cream DAB/SR	14092501	17
Caelo	Wool alcohol cream DAB/SR	14092601	6
Sana	Wool alcohol cream DAB/SR	14092502	1
Euro OTC	Wool alcohol cream DAB/SR	14092601	1
Caelo	Wool alcohol cream DAB/SR	14113201	7
Phönix	Wool alcohol cream DAB/SR	14162502	1
Sanacorp	Wool alcohol cream DAB/SR	14164101	1
Caelo	Wool alcohol cream DAB/SR	14164001	2
Caelo	Wool alcohol cream DAB/SR	14163901	1
Caelo	Wool alcohol cream DAB/SR	14164101	2
Phönix	Wool alcohol cream DAB/SR	14164201	2
Fagron	Wool alcohol cream DAB/SR	14164201	1
Caelo	Wool alcohol cream DAB/SR	14166901ok	1
Phönix	Wool alcohol cream DAB/SR	14169101	3
Caelo	Wool alcohol cream DAB/SR	14169201	4
Phönix	Wool alcohol cream DAB/SR	14170801	1
Caelo	Wool alcohol cream DAB/SR	14170901	1
Euro OTC	Wool alcohol cream DAB/SR	14170901	1
Leyh Pharma/Sanacorp	Wool alcohol cream DAB/SR	142031	1
Caelo	Wool alcohol cream DAB/SR	14216601	7
Caelo	Wool alcohol cream DAB/SR	14216701	3
Caelo	Wool alcohol cream DAB/SR	14216801	5
Caelo	Wool alcohol cream DAB/SR	14216901	2
Caelo	Wool alcohol cream DAB/SR	14216902	15
Fagron	Wool alcohol cream DAB/SR	14216902	1
Caelo	Wool alcohol cream DAB/SR	14216903	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	14246102	1
Caelo	Wool alcohol cream DAB/SR	14246301	7
Phönix	Wool alcohol cream DAB/SR	14246301	1
Caelo	Wool alcohol cream DAB/SR	14246401	6
Caelo	Wool alcohol cream DAB/SR	14246501	9
Gehe	Wool alcohol cream DAB/SR	14246401	1
Ichthyol	Wool alcohol cream DAB/SR	14246501	1
Caelo	Wool alcohol cream DAB/SR	14259401	6
Caelo	Wool alcohol cream DAB/SR	14259402	1
Caelo	Wool alcohol cream DAB/SR	14259501	1
Caelo	Wool alcohol cream DAB/SR	14295301	2
Alliance10.02.2015	Wool alcohol cream DAB/SR	14295502	2
Caelo	Wool alcohol cream DAB/SR	14295502	5
Noweda	Wool alcohol cream DAB/SR	14295402	1
Caelo	Wool alcohol cream DAB/SR	14295402	1
Euro OTC	Wool alcohol cream DAB/SR	14295402	1
Caelo	Wool alcohol cream DAB/SR	14295701	4
Sanacorp	Wool alcohol cream DAB/SR	14295701	1
Caelo	Wool alcohol cream DAB/SR	14295901	3
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	14295901	1
Caelo	Wool alcohol cream DAB/SR	14296001	7
Caelo	Wool alcohol cream DAB/SR	14356301	4
AHD	Wool alcohol cream DAB/SR	14356301	1
Caelo	Wool alcohol cream DAB/SR	14356302	17
Bombastus	Wool alcohol cream DAB/SR	14356302	5
Caelo	Wool alcohol cream DAB/SR	14356303	1
Caelo	Wool alcohol cream DAB/SR	14356402	4
Caelo	Wool alcohol cream DAB/SR	14356401	1
Caelo	Wool alcohol cream DAB/SR	14356501	2
Caelo	Wool alcohol cream DAB/SR	14356601	4
Caelo	Wool alcohol cream DAB/SR	1436201	1
Caelo	Wool alcohol cream DAB/SR	14362102	3
Alliance Health 06.03.2...	Wool alcohol cream DAB/SR	14362102	1
Phönix	Wool alcohol cream DAB/SR	14362102	1
Caelo	Wool alcohol cream DAB/SR	14362301	6
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	14362301	1
Caelo	Wool alcohol cream DAB/SR	14362401	3
Phoenix Goettin	Wool alcohol cream DAB/SR	14362401	2
Phoenix	Wool alcohol cream DAB/SR	14362501	2
Caelo	Wool alcohol cream DAB/SR	1456302	1
Apomix/Sana	Wool alcohol cream DAB/SR	14A180215	3
Leyh Pharma/Sanacorp	Wool alcohol cream DAB/SR	144951	1
Allianc.Healthc../Apomix	Wool alcohol cream DAB/SR	14A180215	1
Caelo	Wool alcohol cream DAB/SR	15011504	1
Caelo	Wool alcohol cream DAB/SR	15017701	2
Phoenix Goettin	Wool alcohol cream DAB/SR	15017801	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	15017801	1
Caelo	Wool alcohol cream DAB/SR	15022201	3
Leyh Pharma/Sanacorp	Wool alcohol cream DAB/SR	150411	1
Caelo	Wool alcohol cream DAB/SR	15073301	1
Caelo	Wool alcohol cream DAB/SR	15073502	4
Fagron	Wool alcohol cream DAB/SR	15073701	1
Caelo	Wool alcohol cream DAB/SR	151110	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	15073701	1
Caelo	Wool alcohol cream DAB/SR	15132101	10
Caelo	Wool alcohol cream DAB/SR	15132301	3
Caelo	Wool alcohol cream DAB/SR	15149902	2
Caelo	Wool alcohol cream DAB/SR	1515101	1
Caelo	Wool alcohol cream DAB/SR	15174103	7
Caelo	Wool alcohol cream DAB/SR	15174104	16

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	15174105	2
Caelo	Wool alcohol cream DAB/SR	15174201	7
Caelo	Wool alcohol cream DAB/SR	15174301	3
Caelo	Wool alcohol cream DAB/SR	15174302	5
Phönix Gotha	Wool alcohol cream DAB/SR	151831	1
Leyh Pharma GmbH	Wool alcohol cream DAB/SR	151831	1
Leyh Pharma/Sanacorp	Wool alcohol cream DAB/SR	151831	1
Leyh/Noweda	Wool alcohol cream DAB/SR	152321	2
Caelo	Wool alcohol cream DAB/SR	15268/601	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	15266501	1
Caelo	Wool alcohol cream DAB/SR	15268302	13
Anzag	Wool alcohol cream DAB/SR	15268302	1
Caelo	Wool alcohol cream DAB/SR	15268303	2
Caelo	Wool alcohol cream DAB/SR	15268401	3
Caelo	Wool alcohol cream DAB/SR	15268402	3
Caelo	Wool alcohol cream DAB/SR	15268501	5
Caelo	Wool alcohol cream DAB/SR	15268601	12
Noweda	Wool alcohol cream DAB/SR	15268601	1
Caelo	Wool alcohol cream DAB/SR	15361701	1
AHD / Leyh	Wool alcohol cream DAB/SR	153711	1
Leyh /Sanacorp	Wool alcohol cream DAB/SR	154421	1
apomix AHD	Wool alcohol cream DAB/SR	15A160216	1
Allianc. Health. / Apomix	Wool alcohol cream DAB/SR	15A160216	1
Apomix/Gehe	Wool alcohol cream DAB/SR	15A230215	1
apomix / Anzag	Wool alcohol cream DAB/SR	15A160216	1
apomix / Noweda	Wool alcohol cream DAB/SR	15A280617	1
Caelo	Wool alcohol cream DAB/SR	16024301	2
Caelo	Wool alcohol cream DAB/SR	16024401	4
Caelo	Wool alcohol cream DAB/SR	16024403	2
Noweda	Wool alcohol cream DAB/SR	16024403	1
Caelo	Wool alcohol cream DAB/SR	16024501	5
Phoenix	Wool alcohol cream DAB/SR	16024501	1
Caelo	Wool alcohol cream DAB/SR	16024502	1
Noweda	Wool alcohol cream DAB/SR	16024601	1
Caelo	Wool alcohol cream DAB/SR	16024601	5
Leyh Pharma/Noweda	Wool alcohol cream DAB/SR	160551	2
Caelo	Wool alcohol cream DAB/SR	161177701	1
Leyh Pharma/Noweda	Wool alcohol cream DAB/SR	161211	1
Caelo	Wool alcohol cream DAB/SR	16120802	1
Leyh Pharma/Noweda	Wool alcohol cream DAB/SR	161531	1
Caelo	Wool alcohol cream DAB/SR	16177401	8
Gehe	Wool alcohol cream DAB/SR	16177401	2
Alte Apotheke	Wool alcohol cream DAB/SR	16177401	1
Caelo	Wool alcohol cream DAB/SR	16177502	6
PHÖNIX	Wool alcohol cream DAB/SR	16177502	1
hochhaus-Apotheke	Wool alcohol cream DAB/SR	16177601	1
Caelo	Wool alcohol cream DAB/SR	16177602	8
Caelo	Wool alcohol cream DAB/SR	16177601	1
Caelo	Wool alcohol cream DAB/SR	16177603	2
Phoenix	Wool alcohol cream DAB/SR	16177701	1
Caelo	Wool alcohol cream DAB/SR	16177701	3
Caelo	Wool alcohol cream DAB/SR	16177801	1
Caelo	Wool alcohol cream DAB/SR	16198401	1
Fagron	Wool alcohol cream DAB/SR	16199003	1
Caelo	Wool alcohol cream DAB/SR	16199401	1
Leyh Pharma/Noweda	Wool alcohol cream DAB/SR	162411	1
Leyh Pharma GmbH, 15.08...	Wool alcohol cream DAB/SR	162411	1
Leyh /Sanacorp	Wool alcohol cream DAB/SR	162631	1
Leyh Pharma/Noweda	Wool alcohol cream DAB/SR	162631	2
Allianz Healthcare	Wool alcohol cream DAB/SR	16312901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix	Wool alcohol cream DAB/SR	16341001	1
phönix	Wool alcohol cream DAB/SR	16344101	1
Caelo	Wool alcohol cream DAB/SR	16344102	2
Caelo	Wool alcohol cream DAB/SR	16344201	8
Caelo	Wool alcohol cream DAB/SR	16344202	1
Gehe	Wool alcohol cream DAB/SR	16344202	1
Caelo	Wool alcohol cream DAB/SR	16344301	10
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	16344301	1
Caelo	Wool alcohol cream DAB/SR	16345801	2
Caelo	Wool alcohol cream DAB/SR	17039202	1
Leyh /Sanacorp	Wool alcohol cream DAB/SR	170231	1
Leyh Pharma/Noweda	Wool alcohol cream DAB/SR	170651	1
Caelo	Wool alcohol cream DAB/SR	17094002	3
Caelo	Wool alcohol cream DAB/SR	17094201	5
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	17094201	1
Caelo	Wool alcohol cream DAB/SR	17104702	2
Caelo	Wool alcohol cream DAB/SR	171128ECA/17149002	1
Caelo	Wool alcohol cream DAB/SR	17178301	7
Caelo	Wool alcohol cream DAB/SR	17178201	3
Caelo	Wool alcohol cream DAB/SR	17178402	1
Caelo	Wool alcohol cream DAB/SR	171851	1
Caelo	Wool alcohol cream DAB/SR	17187601	1
Caelo	Wool alcohol cream DAB/SR	17187701	1
Caelo	Wool alcohol cream DAB/SR	17188001	1
Caelo	Wool alcohol cream DAB/SR	17197301	2
Caelo	Wool alcohol cream DAB/SR	17197302	5
Caelo	Wool alcohol cream DAB/SR	17197903	2
Caelo	Wool alcohol cream DAB/SR	17227401	3
Caelo	Wool alcohol cream DAB/SR	17228702	1
leyh Pharma	Wool alcohol cream DAB/SR	172411	1
Leyh Pharma/Sanacorp	Wool alcohol cream DAB/SR	172411	1
Apomix/Noweda	Wool alcohol cream DAB/SR	17A230215	1
Caelo	Wool alcohol cream DAB/SR	1728701	1
Bombastus	Wool alcohol cream DAB/SR	18000762	1
Bombastus	Wool alcohol cream DAB/SR	18000763	1
Bombastus	Wool alcohol cream DAB/SR	18000765	1
Bombastus	Wool alcohol cream DAB/SR	18000764	1
Bombastus	Wool alcohol cream DAB/SR	18000768	1
Bombastus	Wool alcohol cream DAB/SR	18000766	1
Caelo	Wool alcohol cream DAB/SR	18012301	1
18023501	Wool alcohol cream DAB/SR	18023501	1
Caelo	Wool alcohol cream DAB/SR	18069601	1
Caelo	Wool alcohol cream DAB/SR	18086401	1
Sanacorp W.E:26.03.18 ...	Wool alcohol cream DAB/SR	18012501	1
Apomix	Wool alcohol cream DAB/SR	19A010316	1
Caelo	Wool alcohol cream DAB/SR	18094502	1
Pharmachem GmbH Co.KG,2...	Wool alcohol cream DAB/SR	190816	1
Gehe/Apomix	Wool alcohol cream DAB/SR	19A270214	1
apomix	Wool alcohol cream DAB/SR	20A010316	2
Caelo	Wool alcohol cream DAB/SR	20161102-2	1
Apomix / Phönix	Wool alcohol cream DAB/SR	21A280214	2
Apomix/Noweda	Wool alcohol cream DAB/SR	20A270214	1
Sanacorp	Wool alcohol cream DAB/SR	22/AKD2	1
Euro OTC	Wool alcohol cream DAB/SR	221014	1
Caelo	Wool alcohol cream DAB/SR	23101405	1
apomix / Noweda	Wool alcohol cream DAB/SR	23A181017	1
Sanacorp	Wool alcohol cream DAB/SR	25/AKD3	1
Apomix /Noweda	Wool alcohol cream DAB/SR	25a170316	1
Bombastus	Wool alcohol cream DAB/SR	274683	2
Bombastus	Wool alcohol cream DAB/SR	277856	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	279019	2
Bombastus	Wool alcohol cream DAB/SR	279021	1
apomix	Wool alcohol cream DAB/SR	27A180315	1
Bombastus	Wool alcohol cream DAB/SR	28/7175	1
Bombastus	Wool alcohol cream DAB/SR	280988	1
Bombastus	Wool alcohol cream DAB/SR	280391	1
Bombastus	Wool alcohol cream DAB/SR	281503	2
Bombastus	Wool alcohol cream DAB/SR	282748	2
Bombastus	Wool alcohol cream DAB/SR	282749	1
Bombastus	Wool alcohol cream DAB/SR	282751	4
Bombastus	Wool alcohol cream DAB/SR	282753	2
Bombastus	Wool alcohol cream DAB/SR	282754	6
Bombastus	Wool alcohol cream DAB/SR	283574	1
Bombastus	Wool alcohol cream DAB/SR	284579	4
Bombastus	Wool alcohol cream DAB/SR	284577	1
Bombastus	Wool alcohol cream DAB/SR	284580	1
Bombastus	Wool alcohol cream DAB/SR	284450	1
Bombastus	Wool alcohol cream DAB/SR	284636	3
Bombastus	Wool alcohol cream DAB/SR	2846336	1
Caelo	Wool alcohol cream DAB/SR	284638	1
Bombastus	Wool alcohol cream DAB/SR	284638	1
Bombastus	Wool alcohol cream DAB/SR	284943	1
Bombastus	Wool alcohol cream DAB/SR	284951	4
Bombastus	Wool alcohol cream DAB/SR	284947	3
Bombastus	Wool alcohol cream DAB/SR	285109	1
Bombastus	Wool alcohol cream DAB/SR	285817	1
Bombastus	Wool alcohol cream DAB/SR	285867	1
gehe	Wool alcohol cream DAB/SR	285868	1
Bombastus	Wool alcohol cream DAB/SR	285869	1
GEHE/B0mbastus/04.14	Wool alcohol cream DAB/SR	285870	1
Noweda	Wool alcohol cream DAB/SR	285870	1
Bombastus	Wool alcohol cream DAB/SR	285870	1
Bombastus	Wool alcohol cream DAB/SR	285872	2
Bombastus	Wool alcohol cream DAB/SR	285874	2
Gehe	Wool alcohol cream DAB/SR	285873	1
Bombastus	Wool alcohol cream DAB/SR	286233	1
Bombastus	Wool alcohol cream DAB/SR	286234	2
Bombastus	Wool alcohol cream DAB/SR	286235	1
Bombastus	Wool alcohol cream DAB/SR	286236	2
Bombastus	Wool alcohol cream DAB/SR	286237	2
Bombastus	Wool alcohol cream DAB/SR	286238	1
Bombastus	Wool alcohol cream DAB/SR	287085	1
Bombastus	Wool alcohol cream DAB/SR	287093	3
Noweda	Wool alcohol cream DAB/SR	287093	1
Bombastus	Wool alcohol cream DAB/SR	287176	1
Bombastus	Wool alcohol cream DAB/SR	287179	2
Gehe 21.08.14	Wool alcohol cream DAB/SR	287175	1
Sanacorp	Wool alcohol cream DAB/SR	287183	1
Bombastus	Wool alcohol cream DAB/SR	288259	4
Bombastus	Wool alcohol cream DAB/SR	288260	2
Sanacorp	Wool alcohol cream DAB/SR	288264	1
Bombastus	Wool alcohol cream DAB/SR	288261	1
Bombastus	Wool alcohol cream DAB/SR	288350	1
Caelo	Wool alcohol cream DAB/SR	288265	1
Bombastus	Wool alcohol cream DAB/SR	288351	1
Bombastus	Wool alcohol cream DAB/SR	288354	2
Bombastus	Wool alcohol cream DAB/SR	288443	6
Bombastus	Wool alcohol cream DAB/SR	288441	1
Bombastus	Wool alcohol cream DAB/SR	288355	2
Bombastus	Wool alcohol cream DAB/SR	288446	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	288446	1
Bombastus	Wool alcohol cream DAB/SR	289087	1
Bombastus	Wool alcohol cream DAB/SR	289088	2
Caelo	Wool alcohol cream DAB/SR	289089	1
Bombastus	Wool alcohol cream DAB/SR	289091	3
Anzag	Wool alcohol cream DAB/SR	289091	1
Bombastus	Wool alcohol cream DAB/SR	289394	1
Bombastus	Wool alcohol cream DAB/SR	289395	2
Bombastus	Wool alcohol cream DAB/SR	289390	1
Bombastus	Wool alcohol cream DAB/SR	289396	1
Bombastus	Wool alcohol cream DAB/SR	289391	1
Bombastus	Wool alcohol cream DAB/SR	289487	2
Bombastus	Wool alcohol cream DAB/SR	289491	1
Bombastus	Wool alcohol cream DAB/SR	289490	1
Bombastus	Wool alcohol cream DAB/SR	289620	2
Bombastus	Wool alcohol cream DAB/SR	289622	2
Bombastus	Wool alcohol cream DAB/SR	289632	1
Bombastus	Wool alcohol cream DAB/SR	289623	3
Caelo	Wool alcohol cream DAB/SR	29/AKD12	1
Bombastus	Wool alcohol cream DAB/SR	291265	1
Bombastus	Wool alcohol cream DAB/SR	291272	1
Bombastus	Wool alcohol cream DAB/SR	291267	1
Bombastus	Wool alcohol cream DAB/SR	291348	2
Bombastus	Wool alcohol cream DAB/SR	291347	1
Bombastus	Wool alcohol cream DAB/SR	291350	1
Bombastus	Wool alcohol cream DAB/SR	291352	3
Bombastus	Wool alcohol cream DAB/SR	291349	1
Noweda	Wool alcohol cream DAB/SR	291431	1
Bombastus	Wool alcohol cream DAB/SR	291433	1
Bombastus	Wool alcohol cream DAB/SR	291432	1
Bombastus	Wool alcohol cream DAB/SR	291476	4
Bombastus	Wool alcohol cream DAB/SR	291478	1
Bombastus	Wool alcohol cream DAB/SR	291479	3
Bombastus	Wool alcohol cream DAB/SR	291477	1
Noweda	Wool alcohol cream DAB/SR	291478	1
Bombastus	Wool alcohol cream DAB/SR	292693	2
Bombastus	Wool alcohol cream DAB/SR	292694	1
Bombastus	Wool alcohol cream DAB/SR	292769	1
Bombastus	Wool alcohol cream DAB/SR	292695	2
Bombastus	Wool alcohol cream DAB/SR	292697	1
Bombastus	Wool alcohol cream DAB/SR	293063	1
Bombastus	Wool alcohol cream DAB/SR	293064	1
Bombastus	Wool alcohol cream DAB/SR	293069	3
Bombastus	Wool alcohol cream DAB/SR	293067	1
Bombastus	Wool alcohol cream DAB/SR	293066	1
Caelo	Wool alcohol cream DAB/SR	293069	1
Bombastus	Wool alcohol cream DAB/SR	294271	2
Bombastus	Wool alcohol cream DAB/SR	294272	3
Bombastus	Wool alcohol cream DAB/SR	294273	1
Bombastus	Wool alcohol cream DAB/SR	294333	1
Bombastus	Wool alcohol cream DAB/SR	294332	1
Noweda	Wool alcohol cream DAB/SR	294333	2
Bombastus	Wool alcohol cream DAB/SR	294527	2
Bombastus	Wool alcohol cream DAB/SR	294529	1
Noweda	Wool alcohol cream DAB/SR	294529	1
Bombastus	Wool alcohol cream DAB/SR	294532	2
Bombastus	Wool alcohol cream DAB/SR	294533	10
Bombastus	Wool alcohol cream DAB/SR	295218	1
Bombastus	Wool alcohol cream DAB/SR	295462	3
Bombastus	Wool alcohol cream DAB/SR	295463	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	295463	1
Bombastus	Wool alcohol cream DAB/SR	295464	1
Bombastus	Wool alcohol cream DAB/SR	295528	1
Bombastus	Wool alcohol cream DAB/SR	295465	1
Bombastus	Wool alcohol cream DAB/SR	295529	6
Bombastus	Wool alcohol cream DAB/SR	295532	1
Noweda	Wool alcohol cream DAB/SR	295530	2
Bombastus	Wool alcohol cream DAB/SR	295536	5
Bombastus	Wool alcohol cream DAB/SR	295534	1
Euro OTC	Wool alcohol cream DAB/SR	295536	1
Bombastus	Wool alcohol cream DAB/SR	296127	2
Caelo	Wool alcohol cream DAB/SR	296128	1
Bombastus	Wool alcohol cream DAB/SR	296128	1
Bombastus	Wool alcohol cream DAB/SR	296131	1
Caelo	Wool alcohol cream DAB/SR	296131	2
Bombastus	Wool alcohol cream DAB/SR	296132	2
Bombastus	Wool alcohol cream DAB/SR	296133	3
Bombastus	Wool alcohol cream DAB/SR	296134	2
Bombastus	Wool alcohol cream DAB/SR	297411	1
Bombastus	Wool alcohol cream DAB/SR	297413	4
ichthyol	Wool alcohol cream DAB/SR	297414	1
Bombastus	Wool alcohol cream DAB/SR	297523	1
Caelo	Wool alcohol cream DAB/SR	297526	1
Bombastus	Wool alcohol cream DAB/SR	297528	3
Bombastus	Wool alcohol cream DAB/SR	297527	2
Bombastus	Wool alcohol cream DAB/SR	297526	1
Bombastus	Wool alcohol cream DAB/SR	297530	2
Bombastus	Wool alcohol cream DAB/SR	297529	4
Bombastus	Wool alcohol cream DAB/SR	297531	6
Bombastus	Wool alcohol cream DAB/SR	297533	1
Caelo	Wool alcohol cream DAB/SR	298675	2
Noweda	Wool alcohol cream DAB/SR	298673	1
Bombastus	Wool alcohol cream DAB/SR	298676	1
Bombastus	Wool alcohol cream DAB/SR	298674	1
Bombastus	Wool alcohol cream DAB/SR	298668	1
Bombastus	Wool alcohol cream DAB/SR	298672	1
Bombastus	Wool alcohol cream DAB/SR	298677	1
Bombastus	Wool alcohol cream DAB/SR	298678	2
Caelo	Wool alcohol cream DAB/SR	298677	1
Bombastus	Wool alcohol cream DAB/SR	298679	2
Bombastus	Wool alcohol cream DAB/SR	298813	1
Bombastus	Wool alcohol cream DAB/SR	298814	1
Bombastus	Wool alcohol cream DAB/SR	298815	2
Bombastus	Wool alcohol cream DAB/SR	298816	4
Bombastus	Wool alcohol cream DAB/SR	299228	2
Bombastus	Wool alcohol cream DAB/SR	299227	1
Bombastus	Wool alcohol cream DAB/SR	299230	4
Apomix/Noweda	Wool alcohol cream DAB/SR	299230	1
Bombastus	Wool alcohol cream DAB/SR	299231	4
Bombastus	Wool alcohol cream DAB/SR	299232	2
Bombastus	Wool alcohol cream DAB/SR	299233	3
Bombastus	Wool alcohol cream DAB/SR	299234	6
Bombastus	Wool alcohol cream DAB/SR	299780	1
Bombastus	Wool alcohol cream DAB/SR	299782	2
Bombastus	Wool alcohol cream DAB/SR	299783	2
Bombastus	Wool alcohol cream DAB/SR	299784	3
Bombastus	Wool alcohol cream DAB/SR	299785	2
Bombastus	Wool alcohol cream DAB/SR	299786	2
Bombastus	Wool alcohol cream DAB/SR	299787	1
Bombastus	Wool alcohol cream DAB/SR	300576	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	300577	6
Bombastus	Wool alcohol cream DAB/SR	300578	3
Bombastus	Wool alcohol cream DAB/SR	300583	2
Bombastus	Wool alcohol cream DAB/SR	300579	2
Phoenix	Wool alcohol cream DAB/SR	300583	1
Bombastus	Wool alcohol cream DAB/SR	300581	2
Bombastus	Wool alcohol cream DAB/SR	300831	1
Bombastus	Wool alcohol cream DAB/SR	300833	1
Bombastus	Wool alcohol cream DAB/SR	300835	2
Bombastus	Wool alcohol cream DAB/SR	300837	2
Bombastus	Wool alcohol cream DAB/SR	300836	1
Bombastus	Wool alcohol cream DAB/SR	301305	1
Bombastus	Wool alcohol cream DAB/SR	301308	2
Bombastus	Wool alcohol cream DAB/SR	301306	1
Bombastus	Wool alcohol cream DAB/SR	301309	1
Apomix/Noweda	Wool alcohol cream DAB/SR	301311	1
Bombastus	Wool alcohol cream DAB/SR	301310	2
Caesar & Loretz/AHD	Wool alcohol cream DAB/SR	301311	1
Bombastus	Wool alcohol cream DAB/SR	301312	1
Bombastus	Wool alcohol cream DAB/SR	301883	2
Phönix	Wool alcohol cream DAB/SR	301312	1
Bombastus	Wool alcohol cream DAB/SR	301885	1
Bombastus	Wool alcohol cream DAB/SR	301887	5
Bombastus	Wool alcohol cream DAB/SR	301889	2
Bombastus	Wool alcohol cream DAB/SR	301888	1
Bombastus	Wool alcohol cream DAB/SR	301890	1
Bombastus	Wool alcohol cream DAB/SR	301891	4
Noweda	Wool alcohol cream DAB/SR	302030	2
Bombastus	Wool alcohol cream DAB/SR	302029	1
Bombastus	Wool alcohol cream DAB/SR	302030	1
Bombastus	Wool alcohol cream DAB/SR	302031	4
Bombastus	Wool alcohol cream DAB/SR	302032	3
Bombastus	Wool alcohol cream DAB/SR	302033	2
Bombastus	Wool alcohol cream DAB/SR	3020333	1
Bombastus	Wool alcohol cream DAB/SR	302966	5
Bombastus	Wool alcohol cream DAB/SR	302967	1
Bombastus	Wool alcohol cream DAB/SR	302969	2
Bombastus	Wool alcohol cream DAB/SR	30297	1
Bombastus	Wool alcohol cream DAB/SR	302970	3
Bombastus	Wool alcohol cream DAB/SR	302968	1
Bombastus	Wool alcohol cream DAB/SR	302972	5
GEHE	Wool alcohol cream DAB/SR	303523	1
Bombastus	Wool alcohol cream DAB/SR	303525	1
Bombastus	Wool alcohol cream DAB/SR	303528	2
Caelo	Wool alcohol cream DAB/SR	302972	1
Caelo	Wool alcohol cream DAB/SR	303524	1
Bombastus	Wool alcohol cream DAB/SR	303529	1
Bombastus	Wool alcohol cream DAB/SR	304428	1
Bombastus	Wool alcohol cream DAB/SR	304429	3
Bombastus	Wool alcohol cream DAB/SR	304431	1
Bombastus	Wool alcohol cream DAB/SR	304435	7
Bombastus	Wool alcohol cream DAB/SR	304432	1
Caelo	Wool alcohol cream DAB/SR	304435	1
Bombastus	Wool alcohol cream DAB/SR	304436	4
Caelo	Wool alcohol cream DAB/SR	304436	1
Bombastus	Wool alcohol cream DAB/SR	304437	2
GEHE	Wool alcohol cream DAB/SR	305253	1
Bombastus	Wool alcohol cream DAB/SR	305253	4
Bombastus	Wool alcohol cream DAB/SR	305255	4
Bombastus	Wool alcohol cream DAB/SR	305257	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	305256	1
Bombastus	Wool alcohol cream DAB/SR	305260	2
Gehe	Wool alcohol cream DAB/SR	305257	1
Bombastus	Wool alcohol cream DAB/SR	306083	1
Bombastus	Wool alcohol cream DAB/SR	306084	1
Bombastus	Wool alcohol cream DAB/SR	306085	1
Caelo	Wool alcohol cream DAB/SR	306086	2
Bombastus	Wool alcohol cream DAB/SR	306087	1
Bombastus	Wool alcohol cream DAB/SR	306089	1
Bombastus	Wool alcohol cream DAB/SR	306090	4
Bombastus	Wool alcohol cream DAB/SR	306211	2
Bombastus	Wool alcohol cream DAB/SR	306212	2
Bombastus	Wool alcohol cream DAB/SR	306214	3
Bombastus	Wool alcohol cream DAB/SR	306215	2
GEHE	Wool alcohol cream DAB/SR	306215	4
Bombastus	Wool alcohol cream DAB/SR	306217	1
Bombastus	Wool alcohol cream DAB/SR	306218	2
Bombastus	Wool alcohol cream DAB/SR	306816	1
Bombastus	Wool alcohol cream DAB/SR	307284	3
Bombastus	Wool alcohol cream DAB/SR	306820	1
Bombastus	Wool alcohol cream DAB/SR	306818	1
Bombastus	Wool alcohol cream DAB/SR	307286	2
Bombastus	Wool alcohol cream DAB/SR	307613	3
Bombastus	Wool alcohol cream DAB/SR	307616	1
Bombastus	Wool alcohol cream DAB/SR	307617	3
Bombastus	Wool alcohol cream DAB/SR	307618	3
Bombastus	Wool alcohol cream DAB/SR	307619	1
Bombastus	Wool alcohol cream DAB/SR	307621	1
Bombastus	Wool alcohol cream DAB/SR	307622	1
Bombastus	Wool alcohol cream DAB/SR	307624	3
Bombastus	Wool alcohol cream DAB/SR	307623	1
Bombastus	Wool alcohol cream DAB/SR	308301	3
GEHE	Wool alcohol cream DAB/SR	308302	1
Bombastus	Wool alcohol cream DAB/SR	308303	2
GEHE	Wool alcohol cream DAB/SR	308303	2
Bombastus	Wool alcohol cream DAB/SR	308304	1
Bombastus	Wool alcohol cream DAB/SR	308306	1
Bombastus	Wool alcohol cream DAB/SR	308307	3
GEHE	Wool alcohol cream DAB/SR	308308	1
Bombastus	Wool alcohol cream DAB/SR	308309	1
Bombastus	Wool alcohol cream DAB/SR	308310	1
Bombastus	Wool alcohol cream DAB/SR	308313	1
Bombastus	Wool alcohol cream DAB/SR	308314	1
Bombastus	Wool alcohol cream DAB/SR	308315	2
Caelo	Wool alcohol cream DAB/SR	308685	1
Bombastus	Wool alcohol cream DAB/SR	308316	2
Bombastus	Wool alcohol cream DAB/SR	308686	2
Bombastus	Wool alcohol cream DAB/SR	308684	1
Bombastus	Wool alcohol cream DAB/SR	308688	1
Bombastus	Wool alcohol cream DAB/SR	308687	1
Bombastus	Wool alcohol cream DAB/SR	309640	1
Bombastus	Wool alcohol cream DAB/SR	310579	1
Bombastus	Wool alcohol cream DAB/SR	310581	1
Bombastus	Wool alcohol cream DAB/SR	310582	2
Bombastus	Wool alcohol cream DAB/SR	310584	1
Bombastus	Wool alcohol cream DAB/SR	310583	1
Bombastus	Wool alcohol cream DAB/SR	311448	3
Bombastus	Wool alcohol cream DAB/SR	311454	1
Bombastus	Wool alcohol cream DAB/SR	311455	1
Bombastus	Wool alcohol cream DAB/SR	311458	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	311524	2
GEHE	Wool alcohol cream DAB/SR	311525	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	32A120416	1
Apomix	Wool alcohol cream DAB/SR	32A120416	1
Apomix / Phönix	Wool alcohol cream DAB/SR	33A010414	1
apomix/Noweda	Wool alcohol cream DAB/SR	31A110416	1
Apomix 001	Wool alcohol cream DAB/SR	32A120416	1
All. Healthc. /Apomix	Wool alcohol cream DAB/SR	33A150415	1
apomix/Holdermann	Wool alcohol cream DAB/SR	34A250516	2
Aallianc. Healthc. / Ap...	Wool alcohol cream DAB/SR	35A260516	1
Apomix	Wool alcohol cream DAB/SR	36A250516	1
apomix	Wool alcohol cream DAB/SR	34A250516	1
Apomix/Noweda	Wool alcohol cream DAB/SR	37A090414	1
Apomix/gehe	Wool alcohol cream DAB/SR	37A220415	2
Caelo	Wool alcohol cream DAB/SR	37A090414	1
Apomix	Wool alcohol cream DAB/SR	39A050716	2
Euro OTC	Wool alcohol cream DAB/SR	39A150514	1
Apomix/Noweda	Wool alcohol cream DAB/SR	38A050716	1
Raben-Apotheke Mildenau	Wool alcohol cream DAB/SR	38A150514	1
apomix/Phönix	Wool alcohol cream DAB/SR	39A050716	1
Alliance Healthcare/Apo...	Wool alcohol cream DAB/SR	39A150514	2
Apomix	Wool alcohol cream DAB/SR	39A150514	1
Bombastus	Wool alcohol cream DAB/SR	40/14	1
PKH GmbH	Wool alcohol cream DAB/SR	40A080816	1
Apomix PKH GmbH	Wool alcohol cream DAB/SR	40A080816	1
Allianc. Health. /Apomix	Wool alcohol cream DAB/SR	40A080816	1
Apomix	Wool alcohol cream DAB/SR	40A080816	1
Apomix Anzag	Wool alcohol cream DAB/SR	41A150514	2
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	41A230415	1
Phönix	Wool alcohol cream DAB/SR		1
apomix/Noweda	Wool alcohol cream DAB/SR	42a290816	1
apomix / Anzag	Wool alcohol cream DAB/SR	42A290816	1
Euro OTC	Wool alcohol cream DAB/SR	43A041016	1
apomix/Phönix	Wool alcohol cream DAB/SR	43A200515	1
Apomix/Noweda	Wool alcohol cream DAB/SR	44A041016	1
apomix/Holdermann	Wool alcohol cream DAB/SR	45A041016	1
Raben-Apotheke Mildenau	Wool alcohol cream DAB/SR	45A260514	3
Aallianc. Healthc. / Ap...	Wool alcohol cream DAB/SR	46A011116	2
apomix 004	Wool alcohol cream DAB/SR	46A011116	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	46A27515	1
Bombastus	Wool alcohol cream DAB/SR	47/14	1
Apomix / Phönix	Wool alcohol cream DAB/SR	48A110614	1
Gehe/20.08.14/Apomix	Wool alcohol cream DAB/SR	48A110614	1
Allianc. Health. /Apomix	Wool alcohol cream DAB/SR	49A031116	2
Apomix	Wool alcohol cream DAB/SR	49A031116	1
Apomix	Wool alcohol cream DAB/SR	50A011216	1
Apomix/Gehe	Wool alcohol cream DAB/SR	50A110614	1
apomix	Wool alcohol cream DAB/SR	50A011216	1
Lamotte/Noweda	Wool alcohol cream DAB/SR	553700	1
Apomix / Phoenix	Wool alcohol cream DAB/SR	57A100615	1
Apomix/Phoenix	Wool alcohol cream DAB/SR	56A300614	1
PKH GmbH,Phönix	Wool alcohol cream DAB/SR	57A100615	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	57A300614	1
apomix	Wool alcohol cream DAB/SR	58A100615	1
Bombastus	Wool alcohol cream DAB/SR	60089704	1
Apomix	Wool alcohol cream DAB/SR	61A080714	2
Apomix/Noweda	Wool alcohol cream DAB/SR	62A250615	1
Gehe/apomix/23.09.14	Wool alcohol cream DAB/SR	62A090714	1
apomix/Holdermann	Wool alcohol cream DAB/SR	62A090714	1
Apomix	Wool alcohol cream DAB/SR	64A250715	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
apomix 004	Wool alcohol cream DAB/SR	65A060715	1
Apomix/Alliance Healthc...	Wool alcohol cream DAB/SR	66A060715	1
Apomix Anzag	Wool alcohol cream DAB/SR	69A080814	1
apomix / Anzag	Wool alcohol cream DAB/SR	64A300714	1
Alliance Healthcare/Apo...	Wool alcohol cream DAB/SR	69A080814	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	69A120815	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	72-A190814	1
Apomix	Wool alcohol cream DAB/SR	74A130815	1
apomix/Phönix	Wool alcohol cream DAB/SR	74A130815	1
apomix/Noweda	Wool alcohol cream DAB/SR	75A130815	1
Apomix/Noweda	Wool alcohol cream DAB/SR	77A260815	1
Taoasis	Wool alcohol cream DAB/SR	81/12	1
apomix/Holdermann	Wool alcohol cream DAB/SR	78A160913	1
GEHE/Apomix/31.01.14	Wool alcohol cream DAB/SR	82A260913	1
Apomix / Anzag	Wool alcohol cream DAB/SR	81A100915	1
Apomix/Noweda	Wool alcohol cream DAB/SR	81A290914	1
Apomix	Wool alcohol cream DAB/SR	83A290914	2
Apomix	Wool alcohol cream DAB/SR	85A300914	1
Apomix / Phönix	Wool alcohol cream DAB/SR	85A300914	1
Apomix/Gehe	Wool alcohol cream DAB/SR	87A061015	1
Gehe/apomix/17.12.2014	Wool alcohol cream DAB/SR	87A141014	1
Apomix/Sana	Wool alcohol cream DAB/SR	89A291015	1
Allianc. Health. /Apomix	Wool alcohol cream DAB/SR	89A291015	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	89A291015	1
apomix/Noweda	Wool alcohol cream DAB/SR	90A221014	1
Apomix/Sana	Wool alcohol cream DAB/SR	90A291015	1
apomix/Holdermann	Wool alcohol cream DAB/SR	93A111115	1
Apomix	Wool alcohol cream DAB/SR	92A281014	1
Apomix	Wool alcohol cream DAB/SR	92A111115	1
apomix	Wool alcohol cream DAB/SR	96A111115	3
Apomix	Wool alcohol cream DAB/SR	96A181114	4
Sanacorp	Wool alcohol cream DAB/SR		1
Bombastus	Wool alcohol cream DAB/SR	256922	1
Bombastus	Wool alcohol cream DAB/SR	265187	1
Bombastus	Wool alcohol cream DAB/SR	265303	2
Bombastus	Wool alcohol cream DAB/SR	267043	1
Bombastus	Wool alcohol cream DAB/SR	267154	2
Gehe	Wool alcohol cream DAB/SR	267155	2
Bombastus	Wool alcohol cream DAB/SR	267799	1
Bombastus	Wool alcohol cream DAB/SR	268902	1
Bombastus	Wool alcohol cream DAB/SR	269924	1
Bombastus	Wool alcohol cream DAB/SR	269927	1
Gehe	Wool alcohol cream DAB/SR	269933	2
Bombastus	Wool alcohol cream DAB/SR	270033	2
Gehe	Wool alcohol cream DAB/SR	270032	1
Bombastus	Wool alcohol cream DAB/SR	270035	2
Bombastus	Wool alcohol cream DAB/SR	270037	2
Gehe	Wool alcohol cream DAB/SR	270037	1
Bombastus	Wool alcohol cream DAB/SR	270038	1
Bombastus	Wool alcohol cream DAB/SR	270629	1
Gehe	Wool alcohol cream DAB/SR	270720	1
Bombastus	Wool alcohol cream DAB/SR	270721	1
Bombastus	Wool alcohol cream DAB/SR	270723	1
Bombastus	Wool alcohol cream DAB/SR	270725	2
Bombastus	Wool alcohol cream DAB/SR	271812	1
Bombastus	Wool alcohol cream DAB/SR	271811	1
Bombastus	Wool alcohol cream DAB/SR	272127	1
Gehe	Wool alcohol cream DAB/SR	272254	1
Bombastus	Wool alcohol cream DAB/SR	272257	2
Bombastus	Wool alcohol cream DAB/SR	272254	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	273634	2
Caelo	Wool alcohol cream DAB/SR	273635	1
Bombastus	Wool alcohol cream DAB/SR	273637	1
Bombastus	Wool alcohol cream DAB/SR	273636	1
Bombastus	Wool alcohol cream DAB/SR	274347	1
Gehe/PKH	Wool alcohol cream DAB/SR	3A160113	1
Noweda/apomix	Wool alcohol cream DAB/SR	5A290113	2
Pharmachem.	Wool alcohol cream DAB/SR	81111	1
Apomix	Wool alcohol cream DAB/SR	103A291012	3
Apomix	Wool alcohol cream DAB/SR	101A191213	1
Gehe	Wool alcohol cream DAB/SR	103A291012	1
Apomix	Wool alcohol cream DAB/SR	101A291012	1
Caelo	Wool alcohol cream DAB/SR	9011407	1
Apomix	Wool alcohol cream DAB/SR	104A291012	1
apomix/GEHE/04.04.2013	Wool alcohol cream DAB/SR	108A191112	1
Apomix	Wool alcohol cream DAB/SR	109A191112	1
Caelo	Wool alcohol cream DAB/SR	1113E-01371	1
Caelo	Wool alcohol cream DAB/SR	1117601	1
Pharmachem GmbH / Anzag	Wool alcohol cream DAB/SR	110913	1
Caelo	Wool alcohol cream DAB/SR	11226301	1
Caelo	Wool alcohol cream DAB/SR	11226201	1
Caelo	Wool alcohol cream DAB/SR	11226401	2
Phönix Gotha	Wool alcohol cream DAB/SR	11226301	1
Anzag	Wool alcohol cream DAB/SR	11226402	1
Caesar & Loretz GmbH, 4...	Wool alcohol cream DAB/SR	11276102	2
Gehe	Wool alcohol cream DAB/SR	11276102	2
Caelo	Wool alcohol cream DAB/SR	11276102	1
Caelo	Wool alcohol cream DAB/SR	115A171212	1
Caesar & Loretz GmbH, 4...	Wool alcohol cream DAB/SR	12000201	1
Caelo	Wool alcohol cream DAB/SR	12036301	2
Caelo	Wool alcohol cream DAB/SR	12036401	2
Caelo	Wool alcohol cream DAB/SR	12036402	3
Phönix Gotha	Wool alcohol cream DAB/SR	12036403	1
Caelo	Wool alcohol cream DAB/SR	12036403	1
Gehe	Wool alcohol cream DAB/SR	12036501	1
Caesar & Loretz GmbH, 4...	Wool alcohol cream DAB/SR	12036501	1
anzag	Wool alcohol cream DAB/SR	12036601	1
Caelo	Wool alcohol cream DAB/SR	12056402	1
Caelo	Wool alcohol cream DAB/SR	12101602	1
Caelo	Wool alcohol cream DAB/SR	12121901	1
Caesar & Loretz GmbH, 4...	Wool alcohol cream DAB/SR	12123501	1
Caelo	Wool alcohol cream DAB/SR	12240502	1
Anzag, 19.11.2012	Wool alcohol cream DAB/SR	12205601	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	12209101	1
Caelo	Wool alcohol cream DAB/SR	12260401	1
Anzag 28.03.2013	Wool alcohol cream DAB/SR	12260402	1
Caelo	Wool alcohol cream DAB/SR	1226901	1
Anzag 21.11.12	Wool alcohol cream DAB/SR	12269201	1
Caelo	Wool alcohol cream DAB/SR	1230401	1
Noweda	Wool alcohol cream DAB/SR	12307201	1
Krieger	Wool alcohol cream DAB/SR	12307601	1
Noweda	Wool alcohol cream DAB/SR	12308001	1
Sanacorp	Wool alcohol cream DAB/SR	12307701	1
Sanacorp	Wool alcohol cream DAB/SR	12307801	1
Caesar & Loretz GmbH	Wool alcohol cream DAB/SR	12307801	1
Caelo	Wool alcohol cream DAB/SR	12351801	1
Caelo	Wool alcohol cream DAB/SR	12352001	3
Caesar & Loretz GmbH, 4...	Wool alcohol cream DAB/SR	12352001	1
Phönix Gotha	Wool alcohol cream DAB/SR	12352001	1
12370401	Wool alcohol cream DAB/SR	12370401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool alcohol cream DAB/SR	13013202	8
Caelo	Wool alcohol cream DAB/SR	13008203	1
Post-Apotheke	Wool alcohol cream DAB/SR	13013202	2
Caesar&lorenz,Phoenix	Wool alcohol cream DAB/SR	13013202	1
Caelo	Wool alcohol cream DAB/SR	13061801	1
Caesar & Lorenz GmbH	Wool alcohol cream DAB/SR	13061901	1
Phönix Gotha	Wool alcohol cream DAB/SR	13062101	1
Caesar & Lorenz GmbH	Wool alcohol cream DAB/SR	13062201	1
Phönix	Wool alcohol cream DAB/SR	13062301	1
Caelo	Wool alcohol cream DAB/SR	1306804	1
Caelo	Wool alcohol cream DAB/SR	13068501	4
Caelo	Wool alcohol cream DAB/SR	13068503	3
Alliance healthcare	Wool alcohol cream DAB/SR	13068504	6
Caesar&lorenz,Phoenix	Wool alcohol cream DAB/SR	13068504	1
Caelo	Wool alcohol cream DAB/SR	13068504	8
Fagron	Wool alcohol cream DAB/SR	13068601	1
Phoenix	Wool alcohol cream DAB/SR	13068601	1
Caelo	Wool alcohol cream DAB/SR	13081303	1
Gehe	Wool alcohol cream DAB/SR	13145801	1
Euro OTC	Wool alcohol cream DAB/SR	13145801	1
Phönix Gotha	Wool alcohol cream DAB/SR	13172001	1
Caelo	Wool alcohol cream DAB/SR	13172001	8
Caelo	Wool alcohol cream DAB/SR	13172104	1
Caesar&lorenz,Phoenix	Wool alcohol cream DAB/SR	13172103	1
Caelo	Wool alcohol cream DAB/SR	13172105	2
Caelo	Wool alcohol cream DAB/SR	13190202	1
Caelo	Wool alcohol cream DAB/SR	13279601	2
Caelo	Wool alcohol cream DAB/SR	13280601	1
Caelo	Wool alcohol cream DAB/SR	13280801	2
PHOENIX 21.02.2014	Wool alcohol cream DAB/SR	13369402	1
Caesar & Lorenz GmbH	Wool alcohol cream DAB/SR	13373601	1
Pöfneck Pharmachem	Wool alcohol cream DAB/SR	140312	1
Caelo	Wool alcohol cream DAB/SR	14031602	1
Caesar&lorenz,Phoenix	Wool alcohol cream DAB/SR	14092501	2
Caesar&lorenz,Phoenix	Wool alcohol cream DAB/SR	14200702	1
Caelo	Wool alcohol cream DAB/SR	1424501	1
Phoenix	Wool alcohol cream DAB/SR	14A180214	1
P.W.Beyvers GmbH,Phönix	Wool alcohol cream DAB/SR	15A230215	1
Apomix	Wool alcohol cream DAB/SR	19A120213	1
Caelo	Wool alcohol cream DAB/SR	20101112-1	1
Apomix	Wool alcohol cream DAB/SR	20A070212	1
Apomix	Wool alcohol cream DAB/SR	20A120213	1
Caelo	Wool alcohol cream DAB/SR	21105102	1
Caelo	Wool alcohol cream DAB/SR	2112217035	1
Apomix	Wool alcohol cream DAB/SR	21A070212	1
apomix	Wool alcohol cream DAB/SR	22A200213	1
Caelo	Wool alcohol cream DAB/SR	230512	1
Caelo	Wool alcohol cream DAB/SR	23021303	1
Euro OTC	Wool alcohol cream DAB/SR	2507217035	1
Caelo	Wool alcohol cream DAB/SR	25091204	1
Pharmachem Anzag	Wool alcohol cream DAB/SR	260413	2
Pharmachem / Phoenix	Wool alcohol cream DAB/SR	260413	2
Pharmachem/Phönix	Wool alcohol cream DAB/SR	260413	1
Phoenix (Pharmachem Pöf...	Wool alcohol cream DAB/SR	260712	1
PharmaChem	Wool alcohol cream DAB/SR	260712	2
Pharmachem.Pöfneck	Wool alcohol cream DAB/SR	260712	1
Bombastus	Wool alcohol cream DAB/SR	270730	1
Gehe	Wool alcohol cream DAB/SR	272261	1
Bombastus	Wool alcohol cream DAB/SR	274349	1
Bombastus	Wool alcohol cream DAB/SR	273633	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	274628	1
Bombastus	Wool alcohol cream DAB/SR	274629	1
Bombastus	Wool alcohol cream DAB/SR	274631	1
Bombastus	Wool alcohol cream DAB/SR	274630	1
Caelo	Wool alcohol cream DAB/SR	274635	1
Gehe	Wool alcohol cream DAB/SR	274848	1
Bombastus	Wool alcohol cream DAB/SR	274849	1
Gehe	Wool alcohol cream DAB/SR	274684	1
Alliance Healthcare	Wool alcohol cream DAB/SR	275909	1
Gehe	Wool alcohol cream DAB/SR	274686	1
Bombastus	Wool alcohol cream DAB/SR	275913	1
Bombastus	Wool alcohol cream DAB/SR	275917	1
Euro OTC	Wool alcohol cream DAB/SR	275917	1
Bombastus	Wool alcohol cream DAB/SR	275987	1
Bombastus	Wool alcohol cream DAB/SR	275918	1
Caelo	Wool alcohol cream DAB/SR	275988	1
Bombastus	Wool alcohol cream DAB/SR	276212	1
Bombastus	Wool alcohol cream DAB/SR	276216	1
Bombastus	Wool alcohol cream DAB/SR	276219	2
Bombastus	Wool alcohol cream DAB/SR	277093	1
Bombastus	Wool alcohol cream DAB/SR	277095	1
Bombastus	Wool alcohol cream DAB/SR	277097	2
Bombastus	Wool alcohol cream DAB/SR	277096	1
Bombastus	Wool alcohol cream DAB/SR	277204	1
Bombastus	Wool alcohol cream DAB/SR	277770	2
Caelo	Wool alcohol cream DAB/SR	277855	1
Bombastus	Wool alcohol cream DAB/SR	27771	1
Bombastus	Wool alcohol cream DAB/SR	277769	3
Bombastus	Wool alcohol cream DAB/SR	277854	1
Bombastus	Wool alcohol cream DAB/SR	278744	2
Bombastus	Wool alcohol cream DAB/SR	278749	1
Bombastus	Wool alcohol cream DAB/SR	278824	2
Noweda	Wool alcohol cream DAB/SR	278745	1
Bombastus	Wool alcohol cream DAB/SR	278828	2
Bombastus	Wool alcohol cream DAB/SR	278829	1
Noweda	Wool alcohol cream DAB/SR	278829	1
Bombastus	Wool alcohol cream DAB/SR	278831	2
Bombastus	Wool alcohol cream DAB/SR	278832	1
Bombastus	Wool alcohol cream DAB/SR	278833	1
Bombastus	Wool alcohol cream DAB/SR	279014	1
Bombastus	Wool alcohol cream DAB/SR	279020	1
Bombastus	Wool alcohol cream DAB/SR	279870	2
Bombastus	Wool alcohol cream DAB/SR	279874	1
Euro OTC	Wool alcohol cream DAB/SR	279873	1
Bombastus	Wool alcohol cream DAB/SR	279876	1
Bombastus	Wool alcohol cream DAB/SR	279877	2
Gehe	Wool alcohol cream DAB/SR	279876	1
Bombastus	Wool alcohol cream DAB/SR	281504	2
Gehe	Wool alcohol cream DAB/SR	280493	1
Bombastus	Wool alcohol cream DAB/SR	281506	1
Bombastus	Wool alcohol cream DAB/SR	282239	2
Gehe	Wool alcohol cream DAB/SR	282240	1
Gehe	Wool alcohol cream DAB/SR	282242	1
Bombastus	Wool alcohol cream DAB/SR	282242	2
Noweda	Wool alcohol cream DAB/SR	282241	1
Noweda	Wool alcohol cream DAB/SR	282244	1
Bombastus	Wool alcohol cream DAB/SR	282246	2
gehe	Wool alcohol cream DAB/SR	282748	1
Anzag	Wool alcohol cream DAB/SR	282749	1
Bombastus	Wool alcohol cream DAB/SR	283575	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool alcohol cream DAB/SR	284445	1
Bombastus	Wool alcohol cream DAB/SR	284446	1
Bombastus	Wool alcohol cream DAB/SR	284948	2
Gehe/Apomix/25.10.2013	Wool alcohol cream DAB/SR	28A210313	2
Caelo	Wool alcohol cream DAB/SR	29131-01371	1
Apomix	Wool alcohol cream DAB/SR	32A070312	1
Gehe	Wool alcohol cream DAB/SR	32A080413	1
Apomix	Wool alcohol cream DAB/SR	33A070312	1
Apomix	Wool alcohol cream DAB/SR	33A080413	1
Apomix	Wool alcohol cream DAB/SR	39A180413	3
All..Healthc./Apomix	Wool alcohol cream DAB/SR	41A020513	1
apomix	Wool alcohol cream DAB/SR	41A020513	1
Apomix	Wool alcohol cream DAB/SR	41A150515	1
apomix/ANZAG/04.07.2013	Wool alcohol cream DAB/SR	4546233	1
Apomix/Phönix	Wool alcohol cream DAB/SR	47A040613	1
Apomix	Wool alcohol cream DAB/SR	47A040613	1
All..Healthc./Apomix	Wool alcohol cream DAB/SR	47A040613	1
Gehe	Wool alcohol cream DAB/SR	47a040613	1
Apomix	Wool alcohol cream DAB/SR	51A070512	1
Apomix	Wool alcohol cream DAB/SR	56A270613	1
Apomix	Wool alcohol cream DAB/SR	56A300713	1
Apomix	Wool alcohol cream DAB/SR	62A200612	1
Apomix	Wool alcohol cream DAB/SR	69A100712	1
Gehe	Wool alcohol cream DAB/SR	72A270813	1
Apomix/Noweda	Wool alcohol cream DAB/SR	81A260913	1
Apomix	Wool alcohol cream DAB/SR	72A270913	1
apomix PKH GmbH Halle	Wool alcohol cream DAB/SR	88A231013	1
Gehe	Wool alcohol cream DAB/SR	80A010812	1
Apomix	Wool alcohol cream DAB/SR	90A021012	1
Gehe	Wool alcohol cream DAB/SR	95A161012	1
apomix/Anzag	Wool alcohol cream DAB/SR	89A041113	1
Apomix	Wool alcohol cream DAB/SR	96A161012	1
Apomix	Wool alcohol cream DAB/SR	98A241012	1

- 88 135 spectra from 501 *Apo-Ident* customers from a total of 25 959 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Wool alcohol cream DAB/SR* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Wool alcohol cream DAB/SR* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	700	0	34 844
Type B	0	4140	0	89 734
Type C	0	3069	148	88 135

The substance/substance group *Wool alcohol cream DAB/SR* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 99.1429 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 99.8551 %)
Type C	100.0000 % (> 99.6586 %)	95.3994 % (> 95.3062 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31229	32160	6.23	26.81
31248	31287	1.71	28.26
31271	31271	0.00	23.74
31287	31287	0.00	27.02
31447	31447	0.00	28.55
31455	31455	0.00	32.71
31481	31926	2.37	33.42
31923	31923	0.00	28.17
31926	31926	0.00	30.46
31928	31928	0.00	28.24
31930	31930	0.00	22.48
31933	31933	0.00	27.39
31934	31934	0.00	26.29
31937	31937	0.00	22.85
32160	32160	0.00	32.60
33115	33115	0.00	27.78
33140	33140	0.00	28.30

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Wool fat**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30743-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Wool fat; Adeps lanae; Adeps lanae anhydricus

Special notes

When selecting the *Wool fat* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Wool fat*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Wool fat	12122403	30743	40	20120514
Henry Lamotte	Wool fat	5320900	30925	40	20120608
Caelo	Wool fat	14086904	31479	40	20140321
Euro OTC	Wool fat	1403023-01	31533	40	AR-14-FG-012696-01
Caelo	Wool fat	15050902	31742	40	20150223

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Euro OTC	Wool fat	1502024-01	31743	40	AR-16-FG-003762-01
Euro OTC	Wool fat	1601040	31847	40	20160127*
Caelo	Wool fat	15198802	31952	40	20150625
Henry Lamotte	Wool fat	8003158002	33939	40	20170628

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 360 spectra of 9 reference samples from the substance/substance group *Wool fat*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 9 different batches.
- 35 184 spectra from a total of 831 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 340 spectra from 11 reference samples of the substance/substance group *Wool fat*.
- These include spectra of independent samples from 6 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Euro OTC	Wool fat	1403023-01	31533 [†]	20
Caelo	Wool fat	15050902	31742 [†]	20
Euro OTC	Wool fat	1502024-01	31743 [†]	20
Euro OTC	Wool fat	1601040	31847 [†]	20
Caelo	Wool fat	15198802	31952 [†]	20

continued on the next page

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Wool fat	152923	31998	40
Euro OTC	Wool fat	1607001	32459	40
Caelo	Wool fat	16221907	33111	40
Euro OTC	Wool fat	1701051	33245	40
Caelo	Wool fat	173075	34017	40
Euro OTC	Wool fat	1807008-01	34212	40

- 93 534 spectra from a total of 2156 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1854 spectra from 368 *Apo-Ident* customers from 317 batches from the substance/substance group *Wool fat*.
- These include spectra of independent samples from 304 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Wool fat	13306203	12
Caelo	Wool fat	13306205	9
Caelo	Wool fat	14086903	15
VDL 17.04.14;Preis 4,50	Wool fat		1
Caelo	Wool fat	13398202	8
Caelo	Wool fat	14086905	15
Caelo	Wool fat	14086906	18
Caelo	Wool fat	14147901	12
Caelo	Wool fat	14147902	18
AHD;15.12.14;4,82EUR	Wool fat	14238002	2
Caelo	Wool fat	14238004	16
Caelo	Wool fat	14238006	14
Caelo	Wool fat	14304501	22
Caelo	Wool fat	15198802	19
Caelo	Wool fat	15158401	22
Anzag/Beiersdorf	Wool fat	15050902	1
AHD;9.7.15; 4,78EUR	Wool fat	15050902	1
Caelo	Wool fat	15198804	17
Caelo	Wool fat	16003901	21
Caelo	Wool fat	16003905	13
Caelo	Wool fat	16044302	38
Caelo	Wool fat	16163006	15
Caelo	Wool fat	17248103	4
Caelo	Wool fat	16345104	18
Caelo	Wool fat	17307507	1

continued on the next page

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool fat	16221902	6
Caelo	Wool fat	17307502	4
Bombastus	Wool fat	288711	2
Bombastus	Wool fat	281849	1
Henry Lamotte Oils	Wool fat	5787700	1
henry lamotte	Wool fat	5528000	1
henry lamotte	Wool fat	5690800	1
Bomastus direkt	Wool fat	60081303	1
Henry Lamotte Oils GmbH...	Wool fat	8002144001	1
Henry LAmotte Oils GmbH...	Wool fat	8002112002	1
Bombastus	Wool fat	272421	1
Caelo	Wool fat	1206702	1
Herbasinica	Wool fat	5264700	2
henry lamotte oils	Wool fat	5264700	1
Caelo	Wool fat	5375400	2
Caelo	Wool fat	5253400	1
Gehe	Wool fat	12036703	2
Caelo	Wool fat	5264700	3
Caelo	Wool fat	12224405	15
Caelo	Wool fat	12321803	11
Bombastus	Wool fat	260535	1
Caelo	Wool fat	1303403	1
Caelo	Wool fat	13039403	28
Phoenix	Wool fat	13039403	1
Caelo	Wool fat	12321807	2
Caelo	Wool fat	13213106	15
Caelo	Wool fat	13039405	14
Caelo	Wool fat	13213105	9
Caelo	Wool fat	13213104	17
Caelo	Wool fat	13094905	17
Caelo	Wool fat	13306206	18
Caelo	Wool fat	13213107	17
Caelo	Wool fat	13306207	16
Bombastus	Wool fat	281850	3
	Wool fat	5198200	1
Lamotte	Wool fat	5198200	1
Lamotte/Anzag	Wool fat	5264700-3384310110	1
henry lamotte	Wool fat	5375400	2
Ichthyol/Sanacorp	Wool fat	5448200	8
Henry Lamotte Oils	Wool fat	5448200	1
henry lamotte	Wool fat	5520800	1
Caelo	Wool fat	14238008	19
Bombastus	Wool fat	267377	2
Bombastus	Wool fat	254683	1
Bombastus	Wool fat	268668	5
Bombastus	Wool fat	272422	1
Bombastus	Wool fat	273197	3
Caelo	Wool fat	1094641	1
Caelo	Wool fat	4061402	1
Caelo	Wool fat	4032015A	1
Caelo	Wool fat	316E-01056	1
Caelo	Wool fat	2061509	1
Caelo	Wool fat	4101404	1
Caelo	Wool fat	5052014A	1
Caelo	Wool fat	7051401	1
Caelo	Wool fat	612211056	1
Caelo	Wool fat	6102014C	1
Caelo	Wool fat	5081406	1
Caelo	Wool fat	8061304	1
Caelo	Wool fat	6111212	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool fat	8011509	1
Caelo	Wool fat	11016906	1
Caelo	Wool fat	11016905	1
Caelo	Wool fat	11031302	1
Caelo	Wool fat	11097402	1
Caelo	Wool fat	11203902	2
Caelo	Wool fat	11203903	1
Noweda	Wool fat	11203903	2
Cordes	Wool fat	11203908	1
Caelo	Wool fat	11261103	2
Purren Apotheke	Wool fat	11261106	1
Hedinger	Wool fat	11261105	2
Caelo	Wool fat	11261102	1
Caelo	Wool fat	11261106	2
Fagron	Wool fat	11L07-T01	1
Caelo	Wool fat	12006702	1
Caelo	Wool fat	12015603	8
Caelo	Wool fat	12015604	4
Caelo	Wool fat	12015605	5
Caelo	Wool fat	12036702	5
Noweda/ 18.06.12 Caela	Wool fat	12015605	1
Noweda	Wool fat	12036702	1
Phönix 14.08.2012	Wool fat	12036703	1
Caelo	Wool fat	12036703	7
Caesar & Loretz GmbH, 4. . .	Wool fat	12036703	1
Caesar & Loretz GmbH	Wool fat	12036703	1
Caelo	Wool fat	12036704	6
Caelo	Wool fat	1206703	1
Caelo	Wool fat	12074501	2
Caelo	Wool fat	12074502	5
Caelo	Wool fat	12091301	1
Gehe	Wool fat	12122403	1
Caelo	Wool fat	12122403	13
Caesar & Loretz GmbH	Wool fat	12122403	1
Anzag/13.02.2013	Wool fat	12122403	1
Caelo	Wool fat	12122404	18
Sanacorp	Wool fat	12122404	3
Caelo	Wool fat	12122406	12
Caelo	Wool fat	1213E-01056	1
Caelo	Wool fat	12224402	1
Caelo	Wool fat	12224403	6
Caelo	Wool fat	12224404	18
Gehe	Wool fat	12224404	1
vdl	Wool fat	12224404	1
Noweda	Wool fat	12224404	4
Alliance Healthcare	Wool fat	12224404	1
Sanacorp	Wool fat	12224405	1
Kehr	Wool fat	12224405	1
Caelo	Wool fat	1231802	1
Caelo	Wool fat	12321801	2
Caelo	Wool fat	12321802	8
Caelo	Wool fat	12321805	11
Phönix 02.04.2013	Wool fat	12321805	1
cesar&lorenz/Phönix	Wool fat	12321805	1
Caelo	Wool fat	13039402	16
Fagron	Wool fat	12K07-T51-00357	1
Caesar & Loretz GmbH, 4. . .	Wool fat	13039402	1
Fagron	Wool fat	12H31-T01	1
Gehe	Wool fat	13039402	1
gehe 05.09.13	Wool fat	13039402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alliance Healthcare	Wool fat	13039403	1
Gehe	Wool fat	13039403	1
Sanacorp	Wool fat	13039405	1
Caelo	Wool fat	13081206	1
Euro OTC	Wool fat	1305037-01	1
Caelo	Wool fat	13094902	4
Caelo	Wool fat	1309211056	1
Caelo	Wool fat	1309405	1
Caelo	Wool fat	13094903	15
Caelo	Wool fat	13094904	1
Caelo	Wool fat	1321308	1
Caelo	Wool fat	1321305	1
Caelo	Wool fat	1321307	1
Phönix	Wool fat	13213104	1
Sanacorp	Wool fat	13213104	1
Gehe	Wool fat	13213104	1
Caesar & Loretz GmbH	Wool fat	13213104	1
Gehe 11.11.2013	Wool fat	13213104	1
Caesar&Loretz GmbH/Anzag	Wool fat	13213105	1
Noweda	Wool fat	13213107	1
Caelo	Wool fat	13213108	20
Caelo	Wool fat	132131104	1
Caelo	Wool fat	13303208	1
Caelo	Wool fat	13306204	7
gehe	Wool fat	13306204	1
Gehe 26.2.14	Wool fat	13306205	2
Phönix	Wool fat	13306205	1
Caesar & Loretz GmbH	Wool fat	13306205	1
Sanacorp	Wool fat	13306206	1
Caelo	Wool fat	13306208	16
Caelo	Wool fat	1336205	1
Caelo	Wool fat	133062707	1
Caelo	Wool fat	13389202	1
Phönix	Wool fat	1336206	1
Caelo	Wool fat	13398201	15
Phoenix	Wool fat	13398201	1
Euro OTC	Wool fat	13398201	1
Sanacorp	Wool fat	13398202	1
noweda	Wool fat	13398202	1
Caelo	Wool fat	13398291	1
Caelo	Wool fat	14022014A	1
Caelo	Wool fat	140713c	1
Phoenix	Wool fat	13398202	1
Caelo	Wool fat	14086902	1
Euro OTC	Wool fat	1403023-01	1
Caesar & Loretz GmbH	Wool fat	14086903	1
GEHE	Wool fat	14086903	1
Phönix	Wool fat	14086903	1
Caelo	Wool fat	14086904	14
Bergchemie	Wool fat	14086904	1
Sanacorp	Wool fat	14086904	1
Noweda	Wool fat	14086905	1
Berkel AHK	Wool fat	14086905	1
Fagron	Wool fat	14086906	1
Caelo	Wool fat	14086907	8
Caelo	Wool fat	14111402	1
Caelo	Wool fat	14121204	1
Gehe	Wool fat	14147901	1
Euro OTC	Wool fat	14147902	3
Bombastus	Wool fat	14147902	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
12/2015	Wool fat	14147902	1
Caelo	Wool fat	14147903	16
Caelo	Wool fat	14147905	1
Caelo	Wool fat	1423006	1
Caelo	Wool fat	14238001	2
Caelo	Wool fat	14238002	11
Caesar&lorenz,AHC	Wool fat	14238002	1
AHD	Wool fat	14238002	2
Caesar & Loretz GmbH	Wool fat	14238002	1
AHDE 19,05.,2015	Wool fat	14238002	1
Caelo	Wool fat	14238003	6
Sanacorp	Wool fat	14238004	1
Caelo	Wool fat	14238005	6
AHD 29.01.15 EK 7,67	Wool fat	14238006	1
Caelo	Wool fat	14238007	6
Euro OTC	Wool fat	14238007	1
Allianz Healthcare 22.0...	Wool fat	14238008	1
Caelo	Wool fat	14238009	18
Caelo	Wool fat	14238010	7
Caelo	Wool fat	142383	1
caesar	Wool fat	14304501	1
Caesar & Loretz GmbH	Wool fat	14304501	1
Caelo	Wool fat	14304502	12
Fagron	Wool fat	14304502	1
Sanacorp	Wool fat	14304502	1
Caelo	Wool fat	14304504	8
Caelo	Wool fat	1434502	1
Euro OTC	Wool fat	1502024-01	2
Caelo	Wool fat	1505090	1
Fagron	Wool fat	14B14-T03	1
Caelo	Wool fat	15031401	1
Caelo	Wool fat	15041306	1
Caelo	Wool fat	15050902	9
Caesar & Loretz GmbH	Wool fat	15050902	1
Caelo	Wool fat	15050903	15
Caleo	Wool fat	15050903	1
Caesar & Loretz GmbH	Wool fat	15050903	1
Fagron	Wool fat	15050904	1
Caelo	Wool fat	15050904	11
Caesar&Loretz	Wool fat	15050903	1
Caelo	Wool fat	15050906	13
Caelo	Wool fat	151458402	1
Caelo	Wool fat	15158402	39
Caelo	Wool fat	15158405	15
Caelo	Wool fat	15158404	1
Caelo	Wool fat	1518401	1
Caelo	Wool fat	1518405	1
Caesar&Loretz/Noweda	Wool fat	15198802	1
Caesar&Loretz	Wool fat	15198802	1
Noweda/Caesar & Loretz ...	Wool fat	15198802	1
Caelo	Wool fat	15198803	17
EuRho	Wool fat	15198803	2
Caelo	Wool fat	15292309	5
Caelo	Wool fat	15292311	4
Caelo	Wool fat	15292316	10
Phoenix	Wool fat	15292316	1
Fagron	Wool fat	15B27-T27-023866	1
Caelo	Wool fat	15421401	1
Caesar&Loretz GmbH	Wool fat	16003901	1
Noweda	Wool fat	16003901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool fat	16003902	26
Sanocorp WE:18.04.16 EK...	Wool fat	16003902	2
Phönix	Wool fat	16003905	1
Caelo	Wool fat	1603211056	1
Euro OTC	Wool fat	1601040-01	1
Spirig (Galderma)	Wool fat	16044302	1
caesar	Wool fat	16044302	1
Caelo	Wool fat	16044303	21
Caelo	Wool fat	16044304	1
Caelo	Wool fat	16044306	7
Gehe	Wool fat	16044306	1
Henry Lamotte Oils GMBH	Wool fat	160915ALA	1
Caelo	Wool fat	16121308	1
Caelo	Wool fat	16102013D	1
Caelo	Wool fat	1613001	1
Caelo	Wool fat	16163001	12
Caelo	Wool fat	1613006	1
Caesar & Loretz GmbH Hi...	Wool fat	16163001	1
Caelo	Wool fat	16163002	10
Caelo	Wool fat	16163005	1
GEHE	Wool fat	16163006	1
Gehe	Wool fat	16163006	1
GEHE	Wool fat	16221902	1
Caelo	Wool fat	16221903	10
Euro OTC	Wool fat	16221903	2
Caelo	Wool fat	16221904	17
Alte Apotheke	Wool fat	16221904	1
Caelo	Wool fat	16221905	4
GEHE	Wool fat	16221905	2
Cealo/Sanacorp	Wool fat	16221905	2
Caelo	Wool fat	16221906	11
Sanacorp	Wool fat	16221906	1
Gehe	Wool fat	16221906	1
Caelo	Wool fat	16221907	11
Ichthyol-Gesellschaft	Wool fat	16221907	2
Phönix	Wool fat	16221907	1
Phoenix	Wool fat	16221907	1
Phönix/Phönix	Wool fat	16221907	1
Caelo	Wool fat	16267201	21
Sana 31.01.17	Wool fat	16267201	1
Gehe	Wool fat	16267201	1
Caesar&Loretz GmbH,25.0...	Wool fat	16267201	2
Noweda	Wool fat	16267201	1
Caelo	Wool fat	16267204	6
Caelo	Wool fat	16267203	2
Caelo	Wool fat	16328201	1
Caesar & Loretz/Noweda	Wool fat	16345104	1
henry lamotte/gehe	Wool fat	16345104	1
Caelo	Wool fat	16345107	5
Caelo	Wool fat	16345106	2
Caesar & Loretz GmbH	Wool fat	16345106	1
Sanacorp EK:8,10 EUR...	Wool fat	16345107	2
Noweda	Wool fat	16345107	1
Fiebig	Wool fat	16345106	1
Bombastus	Wool fat	17000437	1
Henry Lamotte Oils GMBH...	Wool fat	170303AL/800232001	1
Caelo	Wool fat	1713903	1
Hubertus Apotheke	Wool fat	165221907	1
Caelo	Wool fat	17139301	5
Caelo	Wool fat	17139302	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wool fat	17139303	12
Noweda	Wool fat	17139303	1
Caelo	Wool fat	17167701	2
Caesar und Loretz GmbH	Wool fat	17167701	1
Noweda	Wool fat	17167702	1
Caelo	Wool fat	17167702	2
Caelo	Wool fat	17167704	1
Caelo	Wool fat	17248102	1
chaelo/gehe	Wool fat	17248103	1
Noweda	Wool fat	17248104	2
Fagron	Wool fat	17248106	1
Gehe	Wool fat	17248104	1
Caelo	Wool fat	17248107	5
Caelo	Wool fat	17248104	1
Caelo	Wool fat	17307501	1
Caelo	Wool fat	17307505	4
Caelo	Wool fat	17307503	1
Caelo	Wool fat	17307506	1
Caelo	Wool fat	1739301	1
Henry Lamotte Oils GmbH...	Wool fat	180109AL/8003158003	1
Caelo	Wool fat	1767704	1
Caelo	Wool fat	18021502	1
GEHE	Wool fat	18029902	1
Caelo	Wool fat	18029902	1
Caelo	Wool fat	1914M1056	1
Henry Lamotte Oils GmbH...	Wool fat	180320AL/8003158003	1
Caelo	Wool fat	18041503	1
Caelo	Wool fat	2113U01056	1
Caelo	Wool fat	21071502	1
Caelo	Wool fat	23012015A	1
Caelo	Wool fat	24422111203903	1
Caelo	Wool fat	27091204	1
Caelo	Wool fat	26061205	1
Bombastus	Wool fat	275701	3
Bombastus	Wool fat	275888	2
Bombastus	Wool fat	275889	3
Bombastus	Wool fat	277470	1
Bombastus	Wool fat	279840	2
Bombastus	Wool fat	280167	2
Caelo	Wool fat	2814U-01056	1
Bombastus	Wool fat	283526	3
Caelo	Wool fat	286362	1
Bombastus	Wool fat	286361	1
Bombastus	Wool fat	286362	3
Bombastus	Wool fat	288712	4
Hubertus Apotheke	Wool fat	29022016	1
Caelo	Wool fat	29041407	1
Bombastus	Wool fat	290878	1
Bombastus	Wool fat	290972	1
Henry Lamotte/Anzag	Wool fat	2908900	1
Bombastus	Wool fat	296321	4
Bombastus	Wool fat	294438	2
Bombastus	Wool fat	291339	1
Bombastus	Wool fat	291336	2
Bombastus	Wool fat	295371	1
Caelo	Wool fat	294643	1
Bombastus	Wool fat	297572	1
Bombastus	Wool fat	299568	1
Caelo	Wool fat	30101311	1
Bombastus	Wool fat	301843	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Wool fat	299653	1
Euro OTC	Wool fat	301843	1
Bombastus	Wool fat	301844	2
Bombastus	Wool fat	303926	3
Bombastus	Wool fat	303928	2
Bombastus	Wool fat	303927	1
Bombastus	Wool fat	305553	3
Sanaorp	Wool fat	305553	1
Caelo	Wool fat	4013Q1056	1
Caelo	Wool fat	3513U-01056	1
Caelo	Wool fat	31011302	1
Caelo	Wool fat	3112U-01056	1
Caelo	Wool fat	412062-1-P	1
Caelo	Wool fat	3x13398202	1
Sanacorp	Wool fat	46/AKD1	1
Caelo	Wool fat	46/AKD3	1
Phönix	Wool fat	5100700	1
Klenk	Wool fat	5132100	1
Sanacorp	Wool fat	46/AKD6	1
Noweda	Wool fat	5065900	1
Caelo	Wool fat	5163200	1
Anzag am 08.05.2012	Wool fat	5163200	1
Caelo	Wool fat		2
Lamotte/ Sanacorp	Wool fat	5198200	2
Henry Lamotte Oils	Wool fat	5198200	2
Caelo	Wool fat	5211E-01056	1
Phoenix	Wool fat	5216400	1
Caelo	Wool fat	5212Q-01056	1
Lamotte/Sanacorp	Wool fat	5253400	1
Henry Lamotte Oils	Wool fat	5253400	1
Noweda	Wool fat	5264700	3
Lamotte/Krieger	Wool fat	5264700	2
Henry Lamotte Oils GmbH	Wool fat	5264700	1
Finze	Wool fat	5264700	2
Lamotte/Phönix	Wool fat	5264700	1
phönix, 23.11.2012	Wool fat	5264700	1
Fagron	Wool fat	5264700	1
Lamotte/Otto	Wool fat	5264700	1
Anzag 25.02.13 ek: 14,44EUR	Wool fat	5264700	1
Lamotte Sanacorp	Wool fat	5264700	1
Großhandel	Wool fat	5264700	1
lamotte/anzag	Wool fat	5264700	1
Henry Lamotte	Wool fat	5264700	1
Henry Lamotte Oils/Sana...	Wool fat	5320900	1
Lamotte	Wool fat	5320900	1
Phönix	Wool fat	5375400	1
Gehe/ Lamotte	Wool fat	5375400	1
Noweda	Wool fat	5375400	1
Gehe	Wool fat	5375400	2
Henry Lamotte Oils GmbH	Wool fat	5375400	1
Phoenix	Wool fat	5375400	2
Anzag	Wool fat	5375400	2
Lamotte/Otto	Wool fat	5375400	3
Gehe 02.07.2013	Wool fat	5375400	2
Henry Lamotte/GEHE	Wool fat	5375400	1
Lamotte	Wool fat	5375400	2
Henry Lamotte Oils/ Now...	Wool fat	5375400	1
Noweda	Wool fat	5448200	1
NW/ Lamotte	Wool fat	5375400	1
Henry Lamotte	Wool fat	5448200	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Wool fat	5448200	3
Lamotte Sanacorp	Wool fat	5448200	1
Lamotte/Sanacorp	Wool fat	5448200	1
Lamotte/Phoenix	Wool fat	5448200	1
Henry Lamotte Oils GmbH	Wool fat	5448200	2
Caelo	Wool fat	5448200	1
Großhandel	Wool fat	5448200	1
Henry Lamotte Oils/ Gehe	Wool fat	5520800	2
Lamotte	Wool fat	5520800	1
Henry Lamotte	Wool fat	5520800	1
Lamotte/Nowed	Wool fat	5520800	1
Henry Lamotte Oils/ Now...	Wool fat	5520800	2
Caelo	Wool fat	5520800	4
Henry Lamotte Oils GmBh...	Wool fat	5520800	1
Lamotte/Sanacorp	Wool fat	5520800	1
Phönix	Wool fat	5520800	1
Henry Lamotte Oils	Wool fat	5520800	1
Lamotte/Phoenix	Wool fat	5520800	2
Henry Lamotte/Noweda	Wool fat	5520800	1
Henry Lamotte Oils GmbH	Wool fat	5528000	2
Alliance healthcare	Wool fat	5528000	1
Henry Lamotte Oils GmbH...	Wool fat	5528000	1
Henry Lamotte/Sanacorp	Wool fat	5528000	1
Henry Lamotte Oils GmbH...	Wool fat	5528000	1
Caelo	Wool fat	5528000	1
Lamotte	Wool fat	5528000	3
Henry Lamotte	Wool fat	5528000	1
Lamotte/Jenne 30.10.2014	Wool fat	5528000	1
Lamotte/Fiebig	Wool fat	5528000	1
Lamotte/gehe	Wool fat	5528000	1
GEHE 02.12.2014	Wool fat	5528000	1
Lamotte	Wool fat	5650500	1
Henry Lamotte	Wool fat	5650500	1
Henry Lamotte Oils GmbH...	Wool fat	5650500	1
Henry Lamotte	Wool fat	55908900	1
Henry Lamotte Oils GmbH	Wool fat	5650500	1
Lamotte/Phönix	Wool fat	5650500	1
Henry Lamotte / Anzag	Wool fat	5650500	1
Henry Lamotte / Anzag	Wool fat	5690800	3
Lamotte/Noweda	Wool fat	5690800	2
Henry Lamotte	Wool fat	5690800	2
Caelo	Wool fat	5727500	1
Gehe	Wool fat	5763900	1
Henry Lamotte	Wool fat	5763900	2
Henry Lamotte Oils GmbH	Wool fat	5763900	5
Caelo	Wool fat	5763900	2
Henry Lamotte/Sanacorp	Wool fat	5727500	2
Henry Lamotte Oils GmbH	Wool fat	56908000	1
Henry Lamotte Oils / No...	Wool fat	5763900	1
Lamotte / Gehe	Wool fat	5763900	1
noweda/lamotte	Wool fat	5763900	1
Lamotte/Noweda	Wool fat	5787700	1
Henry Lamotte/Noweda	Wool fat	5787700	1
Lamotte/Otto	Wool fat	5787700	1
Lamotte/Phoenix	Wool fat	5787700	1
Lamotte/Fiebig	Wool fat	5787700	1
Lamotte/Ebert	Wool fat	5787700	1
Henry Lamotte / Phoenix	Wool fat	5860600	2
henry lamotte oils	Wool fat	5860600	2
Lamotte/Fiebig	Wool fat	5860600	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte Oils GmbH	Wool fat	5860600	1
Henry Lamotte	Wool fat	5860600	1
Caelo	Wool fat	5860600	1
HenryLamotte	Wool fat	5860700	1
Lamotte/Sanacorp	Wool fat	5860700	4
Lamotte/Noweda	Wool fat	5860700	1
Henry Lamotte	Wool fat	5908900	2
Henry Lamotte/Kehr	Wool fat	5860700	1
henry lamotte oils Gmh/N	Wool fat	5908900	1
Henry Lamotte/ Gehe	Wool fat	5908900	1
Henry Lamotte/Anzag	Wool fat	5908900	1
Lamotte/Noweda	Wool fat	5909000	2
Henry Lamotte Oils GmbH...	Wool fat	5908900	1
Lamotte	Wool fat	5909000	1
Henry Lamotte	Wool fat	5909000	2
Lamotte/Fiebig	Wool fat	5909000	2
Henry Lamotte/Kehr	Wool fat	5909000	2
Lamotte/Jenne	Wool fat	5909000	1
Henry Lamotte / Phoenix	Wool fat	5909000	1
Lamotte/Gehe	Wool fat	5909000	1
Henry Lamotte	Wool fat	59090000	1
HenryLamotte	Wool fat	5909100	2
Henry Lamotte/Sanacorp	Wool fat	5909100	5
Lamotte	Wool fat	6001861001	1
Lamotte/Anzag	Wool fat	6001861001	2
Henry Lamotte / Noweda	Wool fat	6001861002	1
Henry Lamotte	Wool fat	6001861002	2
Henry Lamotte/Anzag	Wool fat	6001861001	1
Henry Lamotte Oils/Phönix	Wool fat	6001861001	1
Lamotte/Sanacorp	Wool fat	6001861002	2
henry lamotte oils	Wool fat	6001861002	1
Henry Lamotte/Noweda	Wool fat	6001861002	1
Lamotte	Wool fat	6001861002	1
Henry Lamotte / Phoenix	Wool fat	6001861002	1
Hennry Lamotte/Krieger	Wool fat	6001861002	2
Henry Lamotte/Phönix	Wool fat	6001861002	1
Fagron	Wool fat	6001861002	1
lamotte/anzag	Wool fat	6001861003	1
henry lamotte oils	Wool fat	6991861002	1
Henry Lamotte/Sanacorp	Wool fat	60160	1
Henry Lamotte	Wool fat	8002112001	1
Caelo	Wool fat	75798803	1
Henry Lamotte/Anzag	Wool fat	8002112001	3
phönix	Wool fat	8002112001	2
Lamotte	Wool fat	8002112001	1
Caelo	Wool fat	8002112001	1
Fagron	Wool fat	8002112002	1
Henry Lamotte/Noweda	Wool fat	8002112002	2
Lamotte/Fiebig	Wool fat	8002112002	1
Lamotte	Wool fat	8002112002	2
Henry Lamotte / Phoenix	Wool fat	8002112002	1
Lamotte/Gehe	Wool fat	8002112002	1
Henry Lamotte / Noweda	Wool fat	8002112002	1
Lamotte/Noweda	Wool fat	8002112002	1
Henry Lamotte/ Phönix	Wool fat	8002112002	1
Henry Lamotte/VDL	Wool fat	8002112002	1
Henry Lamotte/ Krieger	Wool fat	8002112002	1
Henry Lamotte Oils/ San...	Wool fat	8002112002	1
phönix	Wool fat	8002112002	1
Henry Lamotte	Wool fat	8002144001	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Henry Lamotte/Noweda	Wool fat	8002144001	1
Henry Lamotte	Wool fat	8002112003	2
Henry Lamotte Oils GmbH	Wool fat	8002144001	1
Henry Lamotte Oils GmbH...	Wool fat	8002144001	2
phönix	Wool fat	8002144001	2
Henry Lamotte/Anzag	Wool fat	8002144001	1
Henry Lamotte Oils GmbH...	Wool fat	8002144001	1
Henry Lamotte Oils / Gehe	Wool fat	8002144001	1
Henry Lamotte Oils GmbH...	Wool fat	8002144001	1
Lamotte/Noweda	Wool fat	8002144002	2
Henry Lamotte/ Noweda	Wool fat	8002320010	4
Henry Lamotte/VDL	Wool fat	8002144002	1
Henry Lamotte	Wool fat	8002320010	1
Henry Lamotte/Phönix	Wool fat	8002320010	1
Henry Lamotte Oils	Wool fat	8002320010	1
Lamotte/Noweda	Wool fat	8002320010	2
Henry Lamotte/ Phoenix	Wool fat	8002320010	1
Sanacorp	Wool fat	8002320010	1
Henry Lamotte Oils GmbH...	Wool fat		1
Lamotte/Gehe	Wool fat	8002320011	1
lamotte/ Phoenix	Wool fat	8002320011	1
Henry Lamotte/Anzag	Wool fat	800244001	1
Henry Lamotte	Wool fat	8002760001	2
HenryLamotte	Wool fat	8002320011	1
Henry Lamotte/Alliance	Wool fat	8002760001	1
Henry Lamotte/Noweda	Wool fat	8002760001	2
Henry Lamotte Oils/ Gehe	Wool fat	8002760002	1
Lamotte/Noweda	Wool fat	8002760002	1
Caelo	Wool fat	82325376	1
Henry Lamotte /Phoenix	Wool fat	8003158002	1
Henry Lamotte/Noweda	Wool fat	8003158002	1
Caelo	Wool fat	AKD17	1
Sanacorp	Wool fat		1
ACM	Wool fat	KNC1879J13	2

- 89 498 spectra from 501 *Apo-Ident* customers from a total of 26 698 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Wool fat* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Wool fat* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	360	0	35 184
Type B	0	340	0	93 534
Type C	0	1831	23	89 498

The substance/substance group *Wool fat* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9703 %)	100.0000 % (> 98.3333 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 98.2353 %)
Type C	100.0000 % (> 99.6586 %)	98.7594 % (> 98.5976 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30743	31743	4.89	19.13
30925	31847	5.05	19.16
31479	31743	1.96	27.62
31533	31533	0.00	23.81
31742	31743	3.18	22.20
31743	31743	0.00	24.40
31847	31847	0.00	24.34
31952	31743	2.03	22.90
33939	31743	3.59	22.95

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **Wormwood oil**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30846-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Wormwood oil; Oleum absinthii

Special notes

When selecting the *Wormwood oil* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Wormwood oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Wormwood oil	11346603	30846	40	20111208
Caelo	Wormwood oil	13347703	31198	40	1403146
Caelo	Wormwood oil	13347704	31423	40	20131028
Caelo	Wormwood oil	13347704	31424	40	20131028

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Wormwood oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 384 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 120 spectra from 4 reference samples of the substance/substance group *Wormwood oil*.
- These include spectra of independent samples from 2 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Wormwood oil	13347704	31423 [†]	20
Caelo	Wormwood oil	13347704	31424 [†]	20
Caelo	Wormwood oil	15316909	33508	40
Caelo	Wormwood oil	170389	33904	40

- 93 754 spectra from a total of 2164 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Type C Spectra from the field, which were not recorded under the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 6 spectra from 3 *Apo-Ident* customers from 3 batches from the substance/substance group *Wormwood oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Wormwood oil	12197809	1
Voegele	Wormwood oil	2013109402	3
Caelo	Wormwood oil	2017115182	2

- 91 346 spectra from 501 *Apo-Ident* customers from a total of 27 006 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Wormwood oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Wormwood oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	120	0	93 754
Type C	0	6	0	91 346

The substance/substance group *Wormwood oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.6600 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30846	31198	4.89	80.48
31198	31198	0.00	79.57
31423	31198	2.06	80.90
31424	31198	1.60	80.86

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	Ylang-ylang oil
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	31300-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Ylang-ylang oil; Oleum cananga odorata; Oleum ylang-ylang

Special notes

When selecting the *Ylang-ylang oil* substance/substance group, the following information is displayed to the user:

- Advice for the identification of essential oils : Some essential oils are listed in the class APIs & excipients, liquid/semi-solid (with analysis certificate) as well as in the class APIs & excipients, liquid/semi-solid (other).

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Ylang-ylang oil*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Ylang-ylang oil	15101104	31791	40	9203
Caelo	Ylang-ylang oil	15101104	31860	40	20150330
Caelo	Ylang-ylang oil	15101109	32175	40	20150330
Caelo	Ylang-ylang oil	163148	33079	40	20161213*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *Ylang-ylang oil*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 3 different batches.
- 35 384 spectra from a total of 837 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 220 spectra from 6 reference samples of the substance/substance group *Ylang-ylang oil*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Ylang-ylang oil	12040606	31300	40
Caelo	Ylang-ylang oil	14192905	31733	60
Caelo	Ylang-ylang oil	15101104	31791 [†]	20
Caelo	Ylang-ylang oil	15101104	31860 [†]	20
Caelo	Ylang-ylang oil	15101109	32350	40
Caelo	Ylang-ylang oil	16314804	33606	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 654 spectra from a total of 2162 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 1 spectra from 1 *Apo-Ident* customers from 1 batches from the substance/substance group *Ylang-ylang oil*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	Ylang-ylang oil	23011502	1

- 91 351 spectra from 501 *Apo-Ident* customers from a total of 27 008 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Ylang-ylang oil* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Ylang-ylang oil* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	120	100	93 654
Type C	0	0	1	91 351

The substance/substance group *Ylang-ylang oil* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	54.5455 % (> 53.1818 %)
Type C	100.0000 % (> 99.6675 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to

exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31791	31791	0.00	54.78
31860	31791	2.72	56.24
32175	31791	5.76	58.08
33079	33079	0.00	48.90

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group **1,3-butanediol**
Substance class APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date 13/12/2018
Report number 30851-2018-12-13
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

1,3-butanediol; 1,3-Butandiolium; Butylenglycolum

Special notes

When selecting the *1,3-butanediol* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7 *European Pharmacopoeia 8th Edition, Basic Version 2014* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Anhang F [Zusatz zu den Modellen der zweiten Stufe](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *1,3-butanediol*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	1,3-butanediol	12J02-N07	30851	40	1402682
Fagron	1,3-butanediol	14G16-B01-296698	31594	40	20140804

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 80 spectra of 2 reference samples from the substance/substance group *1,3-butanediol*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 2 different batches.
- 35 464 spectra from a total of 838 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a [†]. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 60 spectra from 2 reference samples of the substance/substance group *1,3-butanediol*.
- These include spectra of independent samples from 1 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	1,3-butanediol	14G16-B01-296698	31594 [†]	20
Fagron	1,3-butanediol	15E15-B02-332109	33788	40

- 93 814 spectra from a total of 2165 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 3 spectra from 3 *Apo-Ident* customers from 3 batches from the substance/substance group *1,3-butanediol*.
- These include spectra of independent samples from 3 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Caelo	1,3-butanediol	13L19-B03-291208	1
Fagron	1,3-butanediol	S6526964	1
Fagron	1,3-butanediol	13A17-N03	1

- 91 349 spectra from 501 *Apo-Ident* customers from a total of 27 006 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *1,3-butanediol* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *1,3-butanediol* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	80	0	35 464
Type B	0	60	0	93 814
Type C	0	3	0	91 349

The substance/substance group *1,3-butanediol* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9706 %)	100.0000 % (> 92.5000 %)
Type B	100.0000 % (> 99.9762 %)	100.0000 % (> 90.0000 %)
Type C	100.0000 % (> 99.6615 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested

reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
30851	30851	0.00	54.96
31594	30851	6.39	51.88

VALIDATION REPORT

IdentModul 1.3-2019-01

Validated substance/substance group	2-ethylhexyl laurate
Substance class	APIs & excipients, liquid/semi-solid (with analysis certificate)
Report date	13/12/2018
Report number	30849-2018-12-13
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

2-ethylhexyl laurate; 2-Ethylhexylislauromyristas; Cera liquida

Special notes

When selecting the *2-ethylhexyl laurate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-5416-7	<i>European Pharmacopoeia 8th Edition, Basic Version 2014</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Anhang F	Zusatz zu den Modellen der zweiten Stufe

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *2-ethylhexyl laurate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	2-ethylhexyl lau...	14129303	31373	40	20140429
Caelo	2-ethylhexyl lau...	14129304	31507	40	20140429
Caelo	2-ethylhexyl lau...	153825	31878	40	20151113*
Caelo	2-ethylhexyl lau...	163157	33086	40	20161124*

Validation samples

A total of 220 770 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 160 spectra of 4 reference samples from the substance/substance group *2-ethylhexyl laurate*. These samples are listed above in the section [calibration samples](#). The reference samples originate from 4 different batches.
- 35 384 spectra from a total of 836 batches from further 211 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database.

Samples from batches, from which no spectra have been used for the generation of the database are considered as independent samples. The number of batches from which independent samples have delivered *Type B* spectra for the validation is given below and thus represents the number of independent *Type B* samples. Samples from which a number of the spectra has been used for the generation of the database and other spectra have been incorporated into the validation, are marked with a †. The following applies to the remaining unmarked samples: There was at least one other sample from the same batch (other sales containers, other sample ID), of which reference spectra (*Type A*) were used for the generation of the database.

- 212 spectra from 7 reference samples of the substance/substance group *2-ethylhexyl laurate*.
- These include spectra of independent samples from 5 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Sample ID	Spectra
Pharmachem	2-ethylhexyl laurate	010610	30325	12
apomix	2-ethylhexyl laurate	8013080312	30849	40
Caelo	2-ethylhexyl laurate	14129304	31507 [†]	20
Caelo	2-ethylhexyl laurate	153825	31878 [†]	20
Caelo	2-ethylhexyl laurate	163482	33087	40
Caelo	2-ethylhexyl laurate	17065201	33786	40
Caelo	2-ethylhexyl laurate	180194	33885	40

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

- 93 662 spectra from a total of 2160 batches from further 436 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra from the field, which were not recorded under of the control of *Hiperscan GmbH* and have not been used for the generation of the database. *Hiperscan GmbH* accepts the details provided by the customer regarding the manufacturer and batch number mainly unchecked.

- 179 spectra from 86 *Apo-Ident* customers from 51 batches from the substance/substance group *2-ethylhexyl laurate*.
- These include spectra of independent samples from 49 batches, of which no spectra have been used for the generation of the database.

Supplier	Substance	Batch	Spectra
Pharmachem/ Noweda	2-ethylhexyl laurate	10515	4
Pharmachem GmbH&Co.KG	2-ethylhexyl laurate	10515	1
Pharmachem/ Sanacorp	2-ethylhexyl laurate	10515	1
Allianc. Health. / Phar...	2-ethylhexyl laurate	10515	1
Pharma Pößneck	2-ethylhexyl laurate	10515	1
Pharmachem	2-ethylhexyl laurate	10515	3
Fagron	2-ethylhexyl laurate	50414	4
Pharmachem/ Noweda	2-ethylhexyl laurate	10616	2
Pharmachem	2-ethylhexyl laurate	50414	5
pharmachem gmbh	2-ethylhexyl laurate	50414	1
Pharmachem/ Noweda	2-ethylhexyl laurate	50414	2
Pharmachem/Sanacorp	2-ethylhexyl laurate	50414	1
Pharmachem./Allianc.H.	2-ethylhexyl laurate	50414	1
Sanacorp/Pharmachem	2-ethylhexyl laurate	50414	2
Caelo	2-ethylhexyl laurate	13252003	1
Caelo	2-ethylhexyl laurate	14182802	1
apomix	2-ethylhexyl laurate	148-8013010616	3
Caelo	2-ethylhexyl laurate	1482802	1
Caelo	2-ethylhexyl laurate	15382601	1
Caelo	2-ethylhexyl laurate	16143802	1
Caelo	2-ethylhexyl laurate	16143801	1
Caelo	2-ethylhexyl laurate	16143803	3
Caelo	2-ethylhexyl laurate	16143804	1
Caelo	2-ethylhexyl laurate	16245001	4
Caelo	2-ethylhexyl laurate	16245004	1
Caelo	2-ethylhexyl laurate	16245003	1
Caelo	2-ethylhexyl laurate	16315701	2
Caelo	2-ethylhexyl laurate	16348201	3
Caelo	2-ethylhexyl laurate	16315703	1
Caelo	2-ethylhexyl laurate	16348202	4
Caesar & Loretz GmbH	2-ethylhexyl laurate	17065201	1
Caelo	2-ethylhexyl laurate	17065201	3
Caeolo / ahd	2-ethylhexyl laurate	17065201	1
Caelo	2-ethylhexyl laurate	17065202	1
Caelo	2-ethylhexyl laurate	170665201	1
EuRho/Noweda	2-ethylhexyl laurate	18019403	1
Caelo	2-ethylhexyl laurate	18121802	1
Apomix/Anzag	2-ethylhexyl laurate	3114U-03641	1
Apomix/Noweda	2-ethylhexyl laurate	8014260413	2
Phoenix/Apomix	2-ethylhexyl laurate	8015010115	1
Apomix/Sanacorp	2-ethylhexyl laurate	8015010115ZA48/15	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	2-ethylhexyl laurate	13018601	1
apomix/Sanacorp	2-ethylhexyl laurate	8015050414	1
Pharmachem GmbH	2-ethylhexyl laurate	131111	2
Caelo	2-ethylhexyl laurate		1
Caelo	2-ethylhexyl laurate	21091303	1
apomix	2-ethylhexyl laurate	8325111111	1
Phönix Gotha	2-ethylhexyl laurate	131111	1
Gehe	2-ethylhexyl laurate	130312	2
Caelo	2-ethylhexyl laurate	8014210112	1
Pharmachem/AHD	2-ethylhexyl laurate	10213	4
Caelo	2-ethylhexyl laurate	130312	1
Caelo	2-ethylhexyl laurate	14129304	1
Pharmachem/Noweda	2-ethylhexyl laurate	10515	5
AHD / Pharmachem Pößneck	2-ethylhexyl laurate	10515	8
Gehe / Pharmachem. Pößn...	2-ethylhexyl laurate	10616	2
Pharmachem/Sanacorp	2-ethylhexyl laurate	10616	2
AHD/Pharmachem	2-ethylhexyl laurate	10616	5
Sanacorp/Pharmachem.	2-ethylhexyl laurate	50414	7
Apomix/Noweda	2-ethylhexyl laurate	148-8013011116	1
Sanacorp/Pharmachem	2-ethylhexyl laurate	290317	2
Apomix / Phönix	2-ethylhexyl laurate	40514	1
Phönix Gotha	2-ethylhexyl laurate	10213	1
Caelo	2-ethylhexyl laurate	8015180213	1
Pharmachem.Pößneck	2-ethylhexyl laurate	10515	3
Pharmachem Pößneck, GEHE	2-ethylhexyl laurate	10515	1
Pharmachem Pößneck/Sana...	2-ethylhexyl laurate	10515	1
Pharmachem	2-ethylhexyl laurate	10616	1
apomix	2-ethylhexyl laurate	2140817	1
Pharmachem Pößneck/Sana...	2-ethylhexyl laurate	50414	1
Pharmachem GmbH	2-ethylhexyl laurate	50414	1
apomix	2-ethylhexyl laurate	148-8013010218	1
apomix/Noweda	2-ethylhexyl laurate	148-8013010616	1
apomix	2-ethylhexyl laurate	1488013010616	1
Phoenix	2-ethylhexyl laurate	148-8013010616	1
PKH Halle,Anzag	2-ethylhexyl laurate	148-8013010616	1
Apomix/Anzag	2-ethylhexyl laurate	148-8013010616	1
Apomix	2-ethylhexyl laurate	148-8013010616	2
Apomix/Phönix	2-ethylhexyl laurate	148-8013010616	1
Apomix/Noweda	2-ethylhexyl laurate	148-8013010616	1
Fagron	2-ethylhexyl laurate	148-8013010616	1
apomix/Noweda	2-ethylhexyl laurate	148-8013011116	1
Apomix	2-ethylhexyl laurate	148-80130106ZA182/16	1
Apomix	2-ethylhexyl laurate	148-8013011116	1
Caelo	2-ethylhexyl laurate	148-8014010815	1
Caelo	2-ethylhexyl laurate	50414	2
Apomix	2-ethylhexyl laurate	8015010115	1
Euro OTC	2-ethylhexyl laurate	50414	2
Apomix/Anzag	2-ethylhexyl laurate	813M-03324	1
apomix/ Phoenix	2-ethylhexyl laurate	8015050414	1
Apomix	2-ethylhexyl laurate	8015050414	1
PKH Halle/Sanacorp	2-ethylhexyl laurate	8013120913	1
apomix	2-ethylhexyl laurate	8015050414	1
Pharmachem/Phönix	2-ethylhexyl laurate	10213	1
Apomix/Noweda	2-ethylhexyl laurate	8015050414ZA96/14	1
Pharmachem GmbH	2-ethylhexyl laurate	10213	2
Pharmachem / Phönix	2-ethylhexyl laurate	10213	1
Pharmachem/ Sanacorp	2-ethylhexyl laurate	130312	1
pharmachem GmbH Pößneck	2-ethylhexyl laurate	10213	1
Pharmachem GmbH	2-ethylhexyl laurate	130312	1
Pharmachem/ Noweda	2-ethylhexyl laurate	130312	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Pharmachem GmbH&Co.KG	2-ethylhexyl laurate	130312	1
Pharmachem Anzag	2-ethylhexyl laurate	130312	2
Apomix	2-ethylhexyl laurate	8014210112	1
Pharmachem	2-ethylhexyl laurate	131111	4
Pharmachem/Noweda	2-ethylhexyl laurate	131111	2
Pharmachem/Gehe	2-ethylhexyl laurate	131111	1
Bombastus	2-ethylhexyl laurate	221010	1

- 91 173 spectra from 501 *Apo-Ident* customers from a total of 26 959 batches from a further 333 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *2-ethylhexyl laurate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *2-ethylhexyl laurate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	160	0	35 384
Type B	0	212	0	93 662
Type C	0	170	9	91 173

The substance/substance group *2-ethylhexyl laurate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9704 %)	100.0000 % (> 96.2500 %)
Type B	100.0000 % (> 99.9761 %)	100.0000 % (> 97.1698 %)
Type C	100.0000 % (> 99.6586 %)	94.9721 % (> 93.2961 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
31373	31878	2.32	9.00
31507	31878	2.41	8.11
31878	31878	0.00	9.09
33086	33086	0.00	6.73

Appendix A: Additional calibration samples (Type A)

Not required.

Appendix B: Additional validation samples (*Type B*)

It is necessary that spectra also enter the validation which cannot be identified with this model. In this manner, it is verified that the model also rejects unknown substances. The spectra for these samples were recorded by *HiperScan GmbH*. They are allocated to the *Type B*. They also include the calibration spectra for other models.

The samples originate from 385 batches. From these 15 887 spectra were recorded. The spectra recorded on independent samples of substances which can be identified with the model are listed respectively in the section *Type B* for the individual substances and do not appear again elsewhere in this list.

Supplier	Substance	Batch	Spectra	Certificate
Anthroposan	Calendula mother...	0131106	40	not necessary
apomix	Benzalkonium chl...	05B150612	40	not necessary
apomix	Betamethasone va...	01B200912	40	not necessary
apomix	Clobatesol propi...	01B200912	40	not necessary
apomix	Dexamethasone 1%...	04A270513	40	not necessary
apomix	Dexamethasone 1%...	07A201011	30	not necessary
apomix	Epinephrine solu...	04A160712	40	not necessary
apomix	Hydrophilic metr...	04B280411	40	from supplier
apomix	Hydrophilic metr...	09B201011	40	from supplier
apomix	Hydrophilic metr...	08B290911	20	not necessary
apomix	Hydroxyethylcell...	12B240611	40	not necessary
apomix	Methylrosanilini...	01A260413	40	not necessary
apomix	Methylrosanilini...	03C200712	80	not necessary
apomix	Methylrosanilini...	05A200712	40	not necessary
apomix	Metronidazole 10...	17A271011	20	not necessary
apomix	Metronidazole 10...	07B110913	40	1403475
apomix	Metronidazole 10...	06B120412	40	from supplier
apomix	Oxytetracycline ...	02B160212	40	from supplier
apomix	Prednisolone 1% ...	05A191011	10	not necessary
apomix	Prednisolone 1% ...	07A031212	40	not necessary
apomix	Prednisolone 1% ...	01A100113	40	not necessary
apomix	Searing paste SR	04B130418	40	not necessary
apomix	Searing paste SR	08B060612	40	not necessary
apomix	Triamcinolone ac...	01A210513	40	not necessary
apomix	Triamcinolone ac...	04A160212	40	from supplier
apomix	Triamcinolone ac...	09C111212	40	not necessary
apomix	Triamcinolone ac...	01B250811	10	not necessary
apomix	Zinc oxide cream...	10B130712	20	from supplier
apomix	Zinc oxide oil S...	10B090712	40	from supplier
apomix	Zinc oxide paste...	01A190218	40	not necessary
apomix	Zinc oxide paste...	31B141117	40	not necessary
Audor	Phenoxyethanol	APDEG4324219XP	40	not necessary
Audor	Phenoxyethanol	APDEG4324219P	60	not necessary
Audor	Phenoxyethanol	APDEG4324219P	71	AR-16-FG-003720-01
Audor	Vitamin F liposo...	AP307176FL	40	not necessary
Audor	Vitamin F liposo...	zu erfassen	40	not necessary
Bene	Tacholiniquin® 1...	350K121	40	from supplier
Bombastus	Bone oil	00238230	50	not necessary
Bombastus	Chinese cinnamon...	263892	40	AR-12-FG-006571-01
Bombastus	Propolis tinctur...	273005	40	not necessary
Bombastus	Propolis tinctur...	273005	40	not necessary
Bombastus	Zinc oxide cream...	00269409	40	not necessary
Bombastus	Zinc oxide cream...	00269349	20	not necessary
Bombastus	Zinc oxide lotio...	00243799	72	not necessary
Bombastus	Zinc oxide oil DAC	00269769	40	from supplier
Bombastus	Zinc oxide oil S...	00269651	40	not necessary
Bombastus	Zinc oxide oil S...	00266696	40	not necessary
Bombastus	Zinc oxide paste...	00256014	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Bombastus	Zinc oxide paste...	00267515	60	from supplier
Bombastus	Zinc paste DAB	00264410	40	from supplier
Caelo	Abwaschbare Salb...	180622	40	20180215*
Caelo	Acorus tincture ...	11355109	40	not necessary
Caelo	Alcoholic liniment	10179425	40	not necessary
Caelo	Alcoholic liniment	13351909	40	not necessary
Caelo	Aloe vera oil	110137-96542	40	not necessary
Caelo	Angelica root ti...	11015908	40	not necessary
Caelo	Arnica oil, fat	11086702	40	from supplier
Caelo	Arnica tincture ...	09358714	44	not necessary
Caelo	Arnica tincture ...	10362805	40	not necessary
Caelo	Aromatic tinctur...	09382716	10	not necessary
Caelo	Aromatic tinctur...	11355009	40	not necessary
Caelo	Basil oil	11035902	40	not necessary
Caelo	Basil oil	11346902	40	not necessary
Caelo	Basil oil	13395803	40	not necessary
Caelo	Basil oil	13153601	40	1404505
Caelo	Belladonna tinct...	10043003	74	not necessary
Caelo	Belladonna tinct...	11322601	40	not necessary
Caelo	Benzalkonium chl...	12033409	40	not necessary
Caelo	Benzalkonium chl...	10190001	40	not necessary
Caelo	Benzoic lard (st...	12063902	40	not necessary
Caelo	Birch tar	10327401	40	1402581
Caelo	Bitter orange ti...	09278014	40	not necessary
Caelo	Bitter orange ti...	11187201	40	not necessary
Caelo	Bitter orange ti...	82474059	40	AR-12-FG-008243-01
Caelo	Bitter tincture ...	11015709	40	not necessary
Caelo	Bitter tincture ...	11380108	60	not necessary
Caelo	Calendula oil, fat	11183103	40	not necessary
Caelo	Calendula oil, fat	14063403	40	not necessary
Caelo	Calendula oil, fat	13287203	40	1403102
Caelo	Calendula tinctu...	10052810	38	not necessary
Caelo	Camphor oil 10%	13317402	60	not necessary
Caelo	Camphor oil 10%	13317402	40	1509155
Caelo	Camphor oil 10%	10203901	40	not necessary
Caelo	Carmelite water	14001202	40	not necessary
Caelo	Carotene (ca. 4m...	12139305	40	1402564
Caelo	Cayenne pepper t...	11102412	40	not necessary
Caelo	Cayenne pepper t...	09165404	30	not necessary
Caelo	Chinchona bark t...	11275507	40	not necessary
Caelo	Cinchona bark ti...	10101502	40	not necessary
Caelo	Cinchona bark ti...	12056901	40	not necessary
Caelo	Cineole	173099	40	20171205*
Caelo	Coal tar solution	09345509	40	not necessary
Caelo	Coal tar solution	90392199	40	F 2009 000 986
Caelo	Coconut oil, har...	10122201	40	not necessary
Caelo	Comfrey root tin...	12054007	80	not necessary
Caelo	Comfrey tincture...	11121406	40	not necessary
Caelo	Dimethyl isosorb...	17169803	40	not necessary
Caelo	Dimethyl isosorb...	15427602	60	not necessary
Caelo	Dimethyl isosorb...	16170904	40	not necessary
Caelo	Dimethyl isosorb...	16170902	12	not required
Caelo	Dimethyl isosorb...	15230402	60	1603485
Caelo	Food colorant b...	12038601	40	not necessary
Caelo	Food colorant g...	10320905	40	not necessary
Caelo	Food colorant g...	12012503	40	not necessary
Caelo	Food colorant r...	11370506	40	not necessary
Caelo	Food colorant y...	12022203	40	not necessary
Caelo	Food colorant y...	12022212	40	not necessary
Caelo	Formic acid 85%	13373301	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Formic acid 85%	11194901	40	1402561
Caelo	Gentian tincture...	09092707	61	not necessary
Caelo	Gentian tincture...	zu erfassen	40	not necessary
Caelo	Gentian tincture...	80932238	40	AR-12-FG-007885-01
Caelo	Ginger oil	163186	40	20161213*
Caelo	Ginger oil	15042113	40	not necessary
Caelo	Ginger tincture ...	11323801	40	not necessary
Caelo	Hawthorn tinctur...	11311816	40	not necessary
Caelo	Horse chestnut s...	11194605	40	1403541
Caelo	Horse chestnut s...	10346909	40	not necessary
Caelo	Hydrophobic tric...	4811E-03565	30	from supplier
Caelo	Hydroxyethylcell...	92968409	20	F 2009 000 984
Caelo	Hydroxyethylcell...	10292702	40	not necessary
Caelo	Hyoscyamus oil	162692	40	20160930*
Caelo	Iodine tincture	09356705	70	not necessary
Caelo	Iron (III) chlor...	11232103	20	not necessary
Caelo	Larch terpentine	161264	60	20160428*
Caelo	Lavender spirit	11232303	40	not necessary
Caelo	Lavender spirit	13300803	42	1403101
Caelo	Lavender spirit	13300803	60	not necessary
Caelo	Lily-of-the-vall...	10322804	40	not necessary
Caelo	Lime water	12055002	40	not necessary
Caelo	Macadamia oil, r...	170556	40	20170302*
Caelo	Macadamia oil, r...	13012102	40	1407349
Caelo	Macadamia oil, r...	14179304	60	not necessary
Caelo	Marian thistle s...	11061106	40	from supplier
Caelo	Matricaria flowe...	10127009	20	not necessary
Caelo	Matricaria flowe...	11269917	40	not necessary
Caelo	Matricaria flowe...	10127005	50	not necessary
Caelo	Matricaria liqui...	12360310	40	not necessary
Caelo	Matricaria liqui...	09220405	41	not necessary
Caelo	Matricaria liqui...	12360306	40	not necessary
Caelo	Matricaria oil, ...	12189803	40	not necessary
Caelo	Matricaria oil, ...	14001103	40	not necessary
Caelo	Matricaria oil, ...	810189803	20	not necessary
Caelo	Melissa leaf tin...	11102807	40	not necessary
Caelo	Melissa leaf tin...	zu erfassen	40	not necessary
Caelo	Melissa leaf tin...	10105412	30	not necessary
Caelo	Myrrh tincture (...)	11260914	30	not necessary
Caelo	Nettle herb liqu...	14176902	40	1501199
Caelo	Nettle root tinc...	11187409	40	not necessary
Caelo	Nettle tincture ...	11108801	40	not necessary
Caelo	Opium tincture, ...	12020302	40	from supplier
Caelo	Opium tincture, ...	13394905	40	0914I-02461
Caelo	Orange flower oi...	10267805	40	not necessary
Caelo	Orange flower oi...	13067707	40	1404598
Caelo	Orange flower oi...	13067705	40	not necessary
Caelo	Orange flower wa...	12238104	40	not necessary
Caelo	Peanut oil, hydr...	13339803	40	not necessary
Caelo	Peanut oil, hydr...	11266907	60	AR-13-FG-010027-01
Caelo	Peppermint spiri...	11185902	70	from supplier
Caelo	Peppermint tinct...	10106502	39	not necessary
Caelo	Peppermint tinct...	11207405	40	not necessary
Caelo	Pepsin wine	12277403	60	not necessary
Caelo	Pepsin wine	15449303	40	not necessary
Caelo	Pepsin wine	15449301	60	not necessary
Caelo	Pepsin wine	14339802	45	1503624
Caelo	Peru balsam, ori...	14223909	60	1703349
Caelo	Peru balsam, ori...	12302401	40	not necessary
Caelo	Phosphoric acid ...	91297199	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Plantain herb li...	11332613	40	not necessary
Caelo	Polyacrylic acid...	173012	40	20171027*
Caelo	Prednisolone ace...	13077701	100	not necessary
Caelo	Primrose root li...	12175307	100	1404596
Caelo	Primrose root li...	13330110	60	not necessary
Caelo	Propolis liquid ...	13347602	40	1404589
Caelo	Propolis liquid ...	11220411	40	not necessary
Caelo	Propolis liquid ...	13347606	60	not necessary
Caelo	Propolis liquid ...	11220406	40	not necessary
Caelo	Propolis tinctur...	11269521	40	not necessary
Caelo	Propolis tinctur...	10305415	21	not necessary
Caelo	Pyroligneous aci...	12023203	40	not necessary
Caelo	Ratanhia tinctur...	zu erfassen	40	not necessary
Caelo	Ratanhia tinctur...	09166415	40	not necessary
Caelo	Rose water	11286002	20	not necessary
Caelo	Rosemary leaf ti...	09135005	40	not necessary
Caelo	Sage leaf liquid...	10305313	40	not necessary
Caelo	Sage leaf liquid...	11372203	40	not necessary
Caelo	Sage tincture (7...	09345013	72	not necessary
Caelo	Salicylic acid c...	153947	40	20151130*
Caelo	Salicylic acid c...	153946	40	20151130*
Caelo	Siam benzoin	395534-70585	40	not necessary
Caelo	Siam benzoin tin...	11239201	60	not necessary
Caelo	Siam benzoin tin...	13169905	42	1405623
Caelo	Spearmint water	11210003	40	not necessary
Caelo	St. John's wort ...	13389218	40	not necessary
Caelo	St. John's wort ...	12076913	40	not necessary
Caelo	St. John's wort ...	11195217	40	not necessary
Caelo	St. John's wort ...	10197305	30	not necessary
Caelo	Standardised cay...	12053902	40	not necessary
Caelo	Standardised cay...	13259103	20	not necessary
Caelo	Texapon® CS paste	14087206	60	1501415
Caelo	Texapon® CS paste	12292001	40	not necessary
Caelo	Texapon® CS paste	16268901	40	not necessary
Caelo	Texapon® CS paste	12292006	40	not necessary
Caelo	Texapon® CS paste	16268901	40	not necessary
Caelo	Theriac	10122603	43	not necessary
Caelo	Thyme liquid ext...	10198601	41	not necessary
Caelo	Thyme liquid ext...	13273613	60	not necessary
Caelo	Thyme tincture (...)	10043113	40	not necessary
Caelo	Tormentil tinctu...	10067310	40	not necessary
Caelo	Tormentil tinctu...	11278301	40	not necessary
Caelo	Triamcinolone ac...	12250701	80	not necessary
Caelo	Triamcinolone ac...	11258801	40	not necessary
Caelo	Triamcinolone ac...	11258801	30	not necessary
Caelo	Valerian root li...	13343703	60	not necessary
Caelo	Valerian tinctur...	12113006	40	not necessary
Caelo	Valerian tinctur...	92588399	40	AR-12-FG-007325-02
Caelo	Valerian tinctur...	11254304	40	not necessary
Caelo	Valerian tinctur...	12196203	40	not necessary
Caelo	Valerian tinctur...	12113006	40	not necessary
Caelo	Vitamin A-palmit...	11178101	40	not necessary
Caelo	Vitamin A-palmit...	11178101	40	not necessary
Caelo	Vitamin A-palmit...	12132910	60	not necessary
Caelo	Witch hazel bark...	09371313	40	not necessary
Caelo	Witch hazel bark...	13357902	40	1407627
Caelo	Witch hazel bark...	13357902	40	not necessary
Caelo	Witch hazel leaf...	12095402	40	1402497
Caelo	Witch hazel leaf...	12165010	40	not necessary
Caelo	Witch hazel leaf...	10043706	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Wormwood tinctur...	80681198	40	AR-12-FG-007926-01
Caelo	Wormwood tinctur...	12237901	40	not necessary
Caelo	Wormwood tinctur...	10010910	31	not necessary
Caelo	Yarrow herb tinc...	10169006	40	not necessary
Caelo	Zinc ointment DAB	09347002	40	not necessary
Caelo	Zinc ointment DAB	10293402	30	not necessary
Caelo	Zinc oxide lotio...	11263502	40	not necessary
Caelo	Zinc oxide lotio...	11281401	40	from supplier
Caelo	Zinc oxide lotio...	11013001	40	from supplier
Caelo	Zinc oxide lotio...	11263503	40	AR-12-FG-008390-01
Caelo	Zinc oxide lotio...	10340301	40	not necessary
Caelo	Zinc oxide lotio...	12065202	40	not necessary
Caelo	Zinc oxide lotio...	11281401	40	from supplier
Caelo	Zinc oxide lotio...	10340201	33	not necessary
Caelo	Zinc oxide oil DAC	12289802	40	not necessary
Caelo	Zinc oxide oil DAC	12160002	40	not necessary
Caelo	Zinc oxide paste...	00264201	40	not necessary
Caelo	Zinc oxide shaki...	11230101	40	not necessary
Caelo	Zinc paste DAB	11321701	30	not necessary
Caelo	Zinc paste DAB	11290801	40	not necessary
Caelo	Zinc paste DAB	00264409	40	not necessary
Caelo	Zinc paste DAB	12147801	40	not necessary
Caelo	Zinc paste, soft...	11006801	40	not necessary
Caelo	Zinc paste, soft...	12002802	40	not necessary
Caelo	Zinc paste, soft...	12378402	80	not necessary
Euro OTC	Capsaicin liquid...	1205038-02	40	1406349
Euro OTC	Capsaicin liquid...	1205038-01	40	not necessary
Euro OTC	Sodium lactate s...	1201037-01	40	1507187
Euro OTC	Sodium lactate s...	1310007-01	40	AR-17-FG-012458-01
Euro OTC	Sodium lactate s...	1310007-01	40	not necessary
Euro OTC	Sodium lactate s...	1201037-01	60	not necessary
Euro OTC	Sodium lactate s...	1310007-01	60	not necessary
Euro OTC	Sodium lactate s...	1310007-01	24	not required
Fagron	Ammonium bitumin...	10F28-N06	20	not necessary
Fagron	Apricot aroma	11J11-N03	40	not necessary
Fagron	Apricot aroma	13D08-N08	30	1407103
Fagron	Carbomer gel 1%	1205D002	40	from supplier
Fagron	Hydrogen peroxid...	174A270412	40	not necessary
Fagron	Hydrogen peroxid...	12A05-M07	40	not necessary
Fagron	Iron (III) subsu...	17J23-B12-342837	40	not necessary
Fagron	Mesulfen	11G12-N08	40	not necessary
Fagron	Methylrosanilini...	02A190511	40	not necessary
Fagron	Peru balsam, ori...	13F13-N04	40	1502255
Fagron	Rose water	12A04-N04	40	not necessary
Fagron	Tretinoin soluti...	12C08-N14	40	not necessary
Fagron	Tretinoin soluti...	11K09-N03	40	not necessary
Fagron	Tretinoin soluti...	11K09-N03	40	not necessary
Fagron	Tretinoin soluti...	11K09-N03	40	not necessary
Fagron	Zinc paste DAB	12E22-T09	40	from supplier
Gatt-Koller	Camphor spirit (...)	18.00514	40	1809577
Hänseler AG	Potash soap Ph.H...	2017080141	40	20170823
Hänseler AG	Potash soap Ph.H...	2018080637	40	20180906
Hänseler AG	Potash soap Ph.H...	2017120020	40	not necessary
Hedinger	Ammonia solution...	027132	10	from supplier
Hedinger	Ammonia solution...	027138	40	not necessary
Hedinger	Hydrogen peroxid...	301199	40	from supplier
Hedinger	Hydrogen peroxid...	0623498	40	from supplier
Hedinger	Sodium hypochlor...	045226	40	not necessary
Hedinger	Sodium hypochlor...	025339	40	not necessary
Hedinger	Sodium hypochlor...	046359	31	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Henry Lamotte	Zinc oxide lotio...	5375600	40	not necessary
Henry Lamotte	Zinc paste, soft...	5374900	30	not necessary
Henry Lamotte	Zinc paste, soft...	5067100	40	from supplier
HerbaSinica	%Glycerol + prop...	002	80	not necessary
HerbaSinica	%Glycerol + prop...	001	80	not necessary
Ichthyol	Betamethason-V 1...	09D034	50	AR-12-FG-003630-01
Ichthyol	Clobetasol 0.5% ...	09D035	50	AR-12-FG-003627-01
Ichthyol	Clobetasol 0.5% ...	11D034	40	not necessary
Ichthyol	Gentamicin 1% Co...	09D002	50	AR-12-FG-003628-01
Ichthyol	Ichthyol [®] raw ma...	11A075	20	not necessary
Ichthyol	Methoxsalen 0.00...	11D033	40	from supplier
Infectopharm	Hydrophobic tric...	1410F-03565	30	from supplier
Infectopharm	Hydrophobic tric...	S071003.1	70	from supplier
Infectopharm	Hydrophobic tric...	S071007.1	40	not necessary
Infectopharm	Hydrophobic tric...	S091107.1	40	not necessary
Infectopharm	Hydrophobic tric...	S091106.1	40	not necessary
Infectopharm	Hydrophobic tric...	S091106.1	40	from supplier
Kreuz-Apotheke	Aristarom Peach	60644	40	not necessary
Leyh Pharma	Methylrosanilini...	131711	40	not necessary
Madaus	Tormentil tinctu...	L0900294	30	not necessary
Maros	Aaron's rod moth...	U3029	40	not necessary
Maros	Abrotani mother ...	RU2276	40	not necessary
Maros	Angelica root ti...	06026712	40	not necessary
Maros	Angelica root ti...	06026712	40	not necessary
Maros	Anise fruit tinc...	6856	40	not necessary
Maros	Arnica flower li...	4965	40	not necessary
Maros	Arnica root moth...	U3606	40	not necessary
Maros	Artichoke leaf l...	1009	40	not necessary
Maros	Artichoke leaf t...	5367	40	not necessary
Maros	Ash bark mother ...	U3264	40	not necessary
Maros	Bearberry leaf t...	5377	40	not necessary
Maros	Bilberry leaf mo...	U3749	40	not necessary
Maros	Bitter candytuft...	5408	40	not necessary
Maros	Bittersweet stal...	U3320	40	not necessary
Maros	Black cohosh mot...	U3668	40	not necessary
Maros	Boldo leaf mothe...	U3502	30	not necessary
Maros	Bryony root moth...	U3269	40	not necessary
Maros	Burnet root moth...	U3600	40	not necessary
Maros	Cape aloe tinctu...	02026012	40	not necessary
Maros	Caraway seed tin...	08026612	40	not necessary
Maros	Cistus tincture ...	11025811	40	not necessary
Maros	Common broom mot...	R4688	40	not necessary
Maros	Common wormwood ...	U2448	40	not necessary
Maros	Coneflower herb...	U3041	40	not necessary
Maros	Devil's claw roo...	5573	30	not necessary
Maros	European golden ...	U3057	40	not necessary
Maros	Fennel fruit tin...	6996	40	not necessary
Maros	Figwort mother t...	U3623	40	not necessary
Maros	Ginkgo leaf moth...	U2834	40	not necessary
Maros	Golden rod herb ...	6737	40	not necessary
Maros	Hawthorn fruit m...	U3224	40	not necessary
Maros	Horehound herb m...	U2929	40	not necessary
Maros	Horse chestnut s...	U2882	40	not necessary
Maros	Juniper tincture...	1239	40	not necessary
Maros	Lovage root moth...	U3389	40	not necessary
Maros	Marian thistle f...	U3432	40	not necessary
Maros	Matricaria flowe...	U3347	40	not necessary
Maros	Matricaria liqui...	6615	20	not necessary
Maros	Melissa leaf liq...	6399	40	not necessary
Maros	Mistletoe mother...	U2932	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Maros	Monk's pepper fr...	U3476	40	not necessary
Maros	Motherwort herb ...	U2785	40	not necessary
Maros	Motherwort herb ...	6997	40	not necessary
Maros	Muira puama liqu...	3548	40	not necessary
Maros	Nettle mother ti...	U3633	40	not necessary
Maros	Pansy mother tin...	zu erfassen	40	not necessary
Maros	Pansy mother tin...	40971	40	from supplier
Maros	Pelargonium root...	01028113	40	not necessary
Maros	Peppermint leaf ...	RU2432	40	not necessary
Maros	Primrose root ti...	03026212	40	not necessary
Maros	Primrose root ti...	03026212	40	not necessary
Maros	Ribwort plantain...	6736	40	not necessary
Maros	Robert herb moth...	U3720	40	not necessary
Maros	Rosemary leaf mo...	U3238	40	not necessary
Maros	St. John's wort ...	9712	40	not necessary
Maros	Standardised liq...	1012	20	not necessary
Maros	Strophanthus see...	01028513	40	not necessary
Maros	Sundew herb moth...	U3354	40	not necessary
Maros	Teasel root moth...	U3685	40	not necessary
Maros	Tetterwort mothe...	U3063	40	not necessary
Maros	Thyme liquid ext...	6869	40	not necessary
Maros	Valerian root li...	6398	40	not necessary
Maros	Walnut bark tinc...	1240	40	not necessary
Maros	Waterhorehound h...	U3388	40	not necessary
Maros	Willow bark tinc...	6865	40	not necessary
Maros	Wood violet moth...	U2222	40	not necessary
Maros	Yam mother tinct...	32210	40	not necessary
Otto Fischar	Hydrogen peroxid...	zu erfassen	40	not necessary
Phylak Sachse...	Malva silvestris...	2002012	30	not necessary
Phylak Sachse...	Rhus toxicodendr...	2908014	30	not necessary
Phylak Sachse...	Ruta graveolens ...	2002012	30	not necessary
Ratiopharm	Corticoid-ratiop...	L30263	40	from supplier
Staufen-Pharm	Cocculus spag. z...	06.10.08	30	not necessary
Staufen-Pharm	Oxalic acid D30 ...	09.05.12	80	not necessary
Taoasis	Macadamia oil, r...	4976601-92844	40	not necessary
pending	Dichloroacetic a...	zu erfassen	40	not necessary
pending	Zinc paste, soft...	zu erfassen	40	not necessary

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Appendix C: Additional validation samples (*Type C*)

Spectra for substances from the field which cannot be identified with this model also enter the validation. In this manner, it is verified that the model also rejects unknown substances. The spectra for these samples were recorded by *Apo-Ident* customers. They belong to *Type C*. The information provided by the customer regarding the manufacturer and batch number is taken over by *HiperScan GmbH* to a large extent unchecked.

The samples originate from 2098 batches. From these, 5408 spectra were recorded. The spectra recorded on independent samples of substances from the field which can be identified with the model are listed respectively in the section *Type C* for the individual substances and do not appear again elsewhere in this list.

Supplier	Substance	Batch	Spectra
AHD	Castor oil, hydrated	278661	1
AHD, 16.05.13	Zinc oxide oil DAC	13069701	24
alliance healthcare	Matricaria oil, fat	13204306	2
All..Healthc./Apomix	Hydrophilic metronida...	7B180713	7
Alliance Healthcare	Hydrophobic triclosan...	S091202.1	6
Alliance Healthcare	Zinc paste, soft DAB	281118	1
Alliance Healthcare	Zinc paste, soft DAB	5374900	6
Alliance Healthcare/A...	Hydrophilic metronida...	11201113	1
Alliance Healthcare/A...	Hydrophilic metronida...	6B120613	5
AMH Niemann	Methylrosanilinium ch...	3A060612	1
AMH Niemann GmbH & Co...	Hydrophilic metronida...	4190413	1
AMH Niemann GmbH & Co...	Methylrosanilinium ch...	1b140213	1
AMH Niemann GmbH & Co...	Zinc oxide paste, sof...	25A310712	2
AMH Niemann GmbH& Co...	Hydroxyethylcellulose...	13B310513	1
AMH Niemann GmbH& Co...	Hydroxyethylcellulose...	7B040313	1
AMH Niemann GmbH&Co.KG	Zinc oxide paste, sof...	21B210612	1
AMH Niemann GmbH,Halle	Hydroxyethylcellulose...	10B200412	2
AMH Niemann GmbH+Co Kg	Hydrophilic metronida...	5B170414	1
AMH Niemann/G	Oxytetracycline ointm...	1B200313	1
AMH Niemann/G	Oxytetracycline ointm...	4C300712	1
AMH Niemann/G	Oxytetracycline ointm...	2C160212	1
AMH Niemann/Noweda	Zinc oxide paste SR90	3A231012	1
AMH Niemann/Noweda	Zinc oxide paste SR90	1A150312	1
AMH/ Sanacorp	Zinc oxide paste, sof...	43B150811	1
AMH-Niemann GmbH & Co...	Hydroxyethylcellulose...	2C270112	1
AMH-Niemann GmbH & Co...	Hydroxyethylcellulose...	6C090312	1
AMH-Niemann GmbH & Co...	Hydroxyethylcellulose...	1B270112	1
amv niemann gmbh&co kg	Hydrophilic metronida...	4B060515	1
Anzag	Hydrophobic triclosan...	S021202.1	5
Anzag	Methylrosanilinium ch...	5A200712	1
Anzag	Methylrosanilinium ch...	1A140213	2
Anzag	Opium tincture, stand...	12006212	19
Anzag	Witch hazel bark water	11062810	1
Anzag	Zinc oxide oil SR 90	15B100912	1
Anzag	Zinc oxide paste, sof...	12121501	1
Anzag	Zinc oxide paste, sof...	12207901	1
Anzag	Zinc paste, soft DAB	12069101	6
Anzag	Zinc paste, soft DAB	14023301	4
Anzag am 06.02.2013	Sage tincture (70%) 1...	12196609	4
Anzag 19.12.2012	Castor oil, hydrated	273882	3
Anzag 26.11.12 ek:3,35	Orange flower water	12171703	1
Anzag 27.12.12 ek:5,19	Camphor oil 10%	12269502	10
Anzag 28.11.12	Zinc paste, soft DAB	12208001	9
Anzag 29.04.13 ek: 3,...	Castor oil, hydrated	13096304	1
Anzag ek:10.21EUR 27.08...	Alcoholic liniment	12103301	1
Anzag/ Fischar	Formic acid 85%		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Anzag/Fischar	Hydrogen peroxide sol...	6313035	2
Anzag/Infectopharm	Hydrophobic triclosan...	S091206.1	10
anzag/maros	Opium tincture, stand...	120061212	1
apomix	Benzalkonium chloride...	4a210515	1
apomix	Clobatesol propionate...	1A20213	1
apomix	Hydrogen peroxide sol...	231A160513	1
apomix	Hydrogen peroxide sol...	6A080116	1
apomix	Hydrogen peroxide sol...	6A070114	1
apomix	Hydrophilic metronida...	10C141013	1
apomix	Hydrophilic metronida...	5B020513	6
apomix	Hydrophilic metronida...	9B091013	2
apomix	Hydrophilic metronida...	11B221112	2
apomix	Hydroxyethylcellulose...	20B230712	1
apomix	Hydroxyethylcellulose...	14B310513	1
apomix	Hydroxyethylcellulose...	17B080713	2
apomix	Metronidazole 10% (no...	20B191211	1
apomix	Metronidazole 10% (no...	5B150513	1
apomix	Metronidazole 10% (no...	48170413	1
apomix	Metronidazole 10% (no...	38110313	1
apomix	Searing paste SR	9B280613	2
apomix	Zinc oxide paste SR90	2A230811	1
apomix	Zinc oxide paste, sof...	34A231012	1
Apomix	Benzalkonium chloride...	3A170413	2
Apomix	Benzalkonium chloride...	6B060912	1
Apomix	Hydrogen peroxide sol...	227A160513	1
Apomix	Hydrophilic metronida...	8B110912	5
Apomix	Hydrophilic metronida...	3B140313	4
Apomix	Hydrophilic metronida...	9B041012	2
Apomix	Hydrophilic metronida...	8418213	2
Apomix	Hydrophilic metronida...	11B201113	4
Apomix	Hydrophilic metronida...	4B270412	1
Apomix	Hydrophilic metronida...	2B220213	5
Apomix	Hydrophilic metronida...	10B251012	3
Apomix	Hydrophilic metronida...	2B220212	2
Apomix	Hydrophilic metronida...	1B190112	1
Apomix	Hydrophilic metronida...	1B030113	4
Apomix	Hydroxyethylcellulose...	16B310513	3
Apomix	Hydroxyethylcellulose...	13B260413	3
Apomix	Hydroxyethylcellulose...	11B260413	1
Apomix	Hydroxyethylcellulose...	9B120413	1
Apomix	Hydroxyethylcellulose...	12B260413	1
Apomix	Hydroxyethylcellulose...	6B090312	1
Apomix	Hydroxyethylcellulose...	5B150213	3
Apomix	Hydroxyethylcellulose...	6D230514	1
Apomix	Hydroxyethylcellulose...	27B081113	1
Apomix	Hydroxyethylcellulose...	26B301112	1
Apomix	Hydroxyethylcellulose...	6B040313	2
Apomix	Hydroxyethylcellulose...	14B080612	1
Apomix	Hydroxyethylcellulose...	15B060712	1
Apomix	Hydroxyethylcellulose...	16B080713	1
Apomix	Hydroxyethylcellulose...	21C060913	2
Apomix	Hydroxyethylcellulose...	22B281011	1
Apomix	Hydroxyethylcellulose...	19B141011	1
Apomix	Metronidazole 10% (no...	1B020113	2
Apomix	Oxytetracycline ointm...	1A160212	1
Apomix	Searing paste SR	15B211013	1
Apomix	Zinc oxide cream 25% ...	13A040912	3
Apomix	Zinc oxide cream 25% ...	7B270513	1
Apomix	Zinc oxide oil DAC	26C190911	1
Apomix	Zinc oxide oil SR 90	5B250313	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apomix	Zinc oxide oil SR 90	6B150413	1
Apomix	Zinc oxide paste, sof...	12B140313	1
Apomix	Zinc oxide paste, sof...	3A090113	1
Apomix	Zinc paste, soft DAB	41B051212	1
Apomix	Zinc paste, soft DAB	13A040912	1
ApoMix	Hydrophilic metronida...	6B310512	3
Apomix / 11	Dexamethasone 1% (non...	4A270513	1
Apomix / 11	Triamcinolone acetoni...	5A290713	2
Apomix / 11	Triamcinolone acetoni...	6A260813	2
Apomix / 11	Triamcinolone acetoni...	9C111212	1
Apomix / 11	Triamcinolone acetoni...	10A051213	1
APOMIX / 11	Prednisolone 1% (non...	3A160813	1
Apomix / 20	Searing paste SR	12B150813	2
Apomix / 20	Triamcinolone acetoni...	8A190913	1
Apomix / 30	Triamcinolone acetoni...	7A170913	2
APOMIX / 30	Triamcinolone acetoni...	3A110613	1
Apomix / Phönix	Hydroxyethylcellulose...	21C090813	1
Apomix / Phönix	Hydroxyethylcellulose...	11B120413	1
Apomix Anzag	Hydrophilic metronida...	7B200712	5
Apomix direkt	Benzalkonium chloride...	4A110713	2
Apomix direkt	Benzalkonium chloride...	2B070313	1
APomix Phoenix	Searing paste SR	7B160513	1
apomix PKH	Hydrogen peroxide sol...	472A160913	1
apomix PKH	Hydroxyethylcellulose...	14B080612	1
apomix PKH	Hydroxyethylcellulose...	2B270112	1
apomix PKH	Hydroxyethylcellulose...	9C200412	1
apomix PKH	Hydroxyethylcellulose...	17B060712	1
apomix, Phönix	Zinc oxide paste, sof...	6A110213	2
apomix, PKH GmbH	Methylrosanilinium ch...	2B260413	1
Apomix/ 11	Triamcinolone acetoni...	4A040713	1
apomix/ Gehe	Methylrosanilinium ch...	2B290312	1
Apomix/ Noweda	Hydroxyethylcellulose...	5B090312	1
apomix/ Phoenix	Dexamethasone 1% (non...	2A220313	1
Apomix/ Phönix	Hydroxyethylcellulose...	19B090813	1
Apomix/ Sanacorp	Hydroxyethylcellulose...	19B080713	1
Apomix/ Sanacorp	Zinc oxide cream 25% ...	13B230813	1
Apomix/Anzag	Hydrogen peroxide sol...	67A300113	1
Apomix/Anzag	Hydroxyethylcellulose...	10C120413	1
Apomix/Anzag	Methylrosanilinium ch...	1B260413	1
Apomix/Gehe	Methylrosanilinium ch...	4B200712	1
Apomix/Noweda	Hydrophilic metronida...	8B290813	4
Apomix/Noweda	Hydroxyethylcellulose...	18B080713	2
Apomix/Noweda	Searing paste SR	11B010813	1
Apomix/Noweda	Searing paste SR	3B180213	2
Apomix/Noweda	Zinc oxide cream 25% ...	4A100413	1
apomix/phoenix	Betamethasone valerat...	13C031213	1
apomix/Phönix	Hydrophilic metronida...	11B20113	1
apomix/Phönix	Zinc oxide cream 25% ...	15B230913	1
Apomix/Phönix	Hydrophilic metronida...	12B111213	1
Apomix/Phönix	Zinc oxide paste, sof...	4A090113	3
Apotheke	Hydrogen peroxide sol...	8648	1
Apotheke	Hydrogen peroxide sol...	8784	1
Apotheke	Hydrogen peroxide sol...	7488	1
Apotheke	Hydrogen peroxide sol...	9393	1
Apotheke	Hydrogen peroxide sol...	8350	1
Apotheke	Hydrogen peroxide sol...	8352	1
Apotheke	Hydrogen peroxide sol...	8523	1
Apotheke	Hydrogen peroxide sol...	9949	1
Apotheke	Hydrogen peroxide sol...	9527	1
Apotheke	Hydrogen peroxide sol...	9467	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Apotheke	Hydrogen peroxide sol...	8195	1
Apotheke	Opium tincture, stand...	10007013	17
Audor Pharma	Phenoxyethanol	APDEG432419P	1
Audor Pharma	Phenoxyethanol	APDEG4324219XP	1
Audor Pharma	Vitamin F liposomes	AP207152FL	2
Audor Pharma	Vitamin F liposomes	AP302038FL	2
Audor Pharma	Vitamin F liposomes	AP307176FL	2
Audor Pharma	Zinc oxide lotion, aq...	13139402	12
Berkel AHK	Camphor oil 10%	15028001	3
bom	Zinc paste, soft DAB	286902	1
Bombastus	Castor oil, hydrated	269565	3
Bombastus	Castor oil, hydrated	264468	1
Bombastus	Castor oil, hydrated	275938	1
Bombastus	Castor oil, hydrated	273390	15
Bombastus	Castor oil, hydrated	271754	2
Bombastus	Castor oil, hydrated	270580	5
Bombastus	Castor oil, hydrated	26956512-13	1
Bombastus	Castor oil, hydrated	266593	1
Bombastus	Castor oil, hydrated	266591	1
Bombastus	Castor oil, hydrated	255179	1
Bombastus	Chinese cinnamon oil	269558	1
Bombastus	Chinese cinnamon oil	262660	1
Bombastus	Chinese cinnamon oil	263892	1
Bombastus	Hydroxyethylcellulose...	13190401	14
Bombastus	Propolis tincture 20%	281900	1
Bombastus	Propolis tincture 20%	276954	1
Bombastus	Propolis tincture 20%	275084	1
Bombastus	Propolis tincture 20%	280337	2
Bombastus	St. John's wort oil	272178	1
Bombastus	St. John's wort oil	293138	1
Bombastus	St. John's wort oil	283368	1
Bombastus	St. John's wort oil	282275	1
Bombastus	St. John's wort oil	279926	1
Bombastus	St. John's wort oil	275674	1
Bombastus	Zinc ointment DAB	282182	1
Bombastus	Zinc ointment DAB	275313	4
Bombastus	Zinc ointment DAB	277335	2
Bombastus	Zinc ointment DAB	30071602	1
Bombastus	Zinc oxide cream 25% ...	279248	1
Bombastus	Zinc oxide cream 25% ...	272491	1
Bombastus	Zinc oxide cream 25% ...	276966	2
Bombastus	Zinc oxide cream 25% ...	274259	1
Bombastus	Zinc oxide cream 25% ...	272490	1
Bombastus	Zinc oxide cream 25% ...	281755	2
Bombastus	Zinc oxide lotion, aq...	274933	1
Bombastus	Zinc oxide lotion, aq...	281537	1
Bombastus	Zinc oxide lotion, aq...	277579	1
Bombastus	Zinc oxide lotion, aq...	273199	2
Bombastus	Zinc oxide lotion, aq...	270608	1
Bombastus	Zinc oxide lotion, aq...	284632	1
Bombastus	Zinc oxide oil DAC	301565	1
Bombastus	Zinc oxide oil DAC	278159	7
Bombastus	Zinc oxide oil DAC	282083	1
Bombastus	Zinc oxide oil DAC	274182	1
Bombastus	Zinc oxide oil DAC	273972	1
Bombastus	Zinc oxide oil DAC	273971	1
Bombastus	Zinc oxide oil DAC	273970	1
Bombastus	Zinc oxide oil DAC	269765	1
Bombastus	Zinc oxide oil DAC	269767	1
Bombastus	Zinc oxide oil DAC	269650	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Zinc oxide oil DAC	267450	1
Bombastus	Zinc oxide oil DAC	264960	1
Bombastus	Zinc oxide oil DAC	2608864	1
Bombastus	Zinc oxide oil DAC	269651	3
Bombastus	Zinc oxide oil DAC	279564	1
Bombastus	Zinc oxide oil DAC	286041	2
Bombastus	Zinc oxide oil DAC	263248	1
Bombastus	Zinc oxide oil SR 90	273973	6
Bombastus	Zinc oxide oil SR 90	6858379	1
Bombastus	Zinc oxide oil SR 90	281661	1
Bombastus	Zinc oxide oil SR 90	269651	1
Bombastus	Zinc oxide oil SR 90	277647	1
Bombastus	Zinc oxide paste SR90	279384	1
Bombastus	Zinc oxide paste SR90	282124	1
Bombastus	Zinc oxide paste SR90	279387	1
Bombastus	Zinc oxide paste SR90	276200	4
Bombastus	Zinc oxide paste SR90	276198	1
Bombastus	Zinc oxide paste SR90	273833	1
Bombastus	Zinc oxide paste SR90	256014	1
Bombastus	Zinc oxide paste SR90	266450	1
Bombastus	Zinc oxide paste, sof...	273326	1
Bombastus	Zinc oxide paste, sof...	273329	4
Bombastus	Zinc oxide paste, sof...	273328	3
Bombastus	Zinc oxide paste, sof...	279586	4
Bombastus	Zinc oxide paste, sof...	281586	1
Bombastus	Zinc oxide paste, sof...	271192	2
Bombastus	Zinc oxide paste, sof...	267515	1
Bombastus	Zinc oxide paste, sof...	271196	2
Bombastus	Zinc oxide paste, sof...	273330	3
Bombastus	Zinc oxide paste, sof...	275031	1
Bombastus	Zinc oxide paste, sof...	277662	1
Bombastus	Zinc paste DAB	290176	1
Bombastus	Zinc paste DAB	278270	1
Bombastus	Zinc paste DAB	287850	1
Bombastus	Zinc paste DAB	284662	1
Bombastus	Zinc paste DAB	278272	3
Bombastus	Zinc paste DAB	282037	1
Bombastus	Zinc paste DAB	280732	1
Bombastus	Zinc paste DAB	276604	1
Bombastus	Zinc paste DAB	278865	1
Bombastus	Zinc paste DAB	278273	3
Bombastus	Zinc paste DAB	278861	1
Bombastus	Zinc paste DAB	278864	3
Bombastus	Zinc paste DAB	278863	1
Bombastus	Zinc paste DAB	291718	2
Bombastus	Zinc paste DAB	280733	2
Bombastus	Zinc paste DAB	271136	1
Bombastus	Zinc paste DAB	271135	2
Bombastus	Zinc paste DAB	273890	1
Bombastus	Zinc paste DAB	276146	1
Bombastus	Zinc paste DAB	273894	1
Bombastus	Zinc paste DAB	271134	1
Bombastus	Zinc paste DAB	271133	1
Bombastus	Zinc paste DAB	27113603-15	1
Bombastus	Zinc paste DAB	276145	1
Bombastus	Zinc paste DAB	278274	3
Bombastus	Zinc paste DAB	273893	2
Bombastus	Zinc paste DAB	276143	1
Bombastus	Zinc paste DAB	62/13	1
Bombastus	Zinc paste, soft DAB	282863	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Zinc paste, soft DAB	266381	1
Bombastus	Zinc paste, soft DAB	282862	1
Bombastus	Zinc paste, soft DAB	281122	4
Bombastus	Zinc paste, soft DAB	281121	2
Bombastus	Zinc paste, soft DAB	281119	1
Bombastus	Zinc paste, soft DAB	281116	2
Bombastus	Zinc paste, soft DAB	278229	4
Bombastus	Zinc paste, soft DAB	278227	2
Bombastus	Zinc paste, soft DAB	275948	1
Bombastus	Zinc paste, soft DAB	274973	1
Bombastus	Zinc paste, soft DAB	271341	2
Bombastus	Zinc paste, soft DAB	271344	3
Bombastus	Zinc paste, soft DAB	280233	1
Bombastus	Zinc paste, soft DAB	267590	1
Bombastus	Zinc paste, soft DAB	278226	2
Bombastus	Zinc paste, soft DAB	280235	1
Bombastus	Zinc paste, soft DAB	282859	1
Bombastus	Zinc paste, soft DAB	282861	1
Bombastus	Zinc paste, soft DAB	293485	1
Bombastus	Zinc paste, soft DAB	278225	1
Bombastus	Zinc paste, soft DAB	278223	1
Bombastus	Zinc paste, soft DAB	276539	3
Bombastus	Zinc paste, soft DAB	276538	2
Bombastus	Zinc paste, soft DAB	276537	1
Bombastus	Zinc paste, soft DAB	276536	1
Bombastus	Zinc paste, soft DAB	271340	4
Bombastus	Zinc paste, soft DAB	270516	2
Bombastus	Zinc paste, soft DAB	270465	3
Bombastus	Zinc paste, soft DAB	270461	1
Bombastus	Zinc paste, soft DAB	265457	1
Bombastus	Zinc paste, soft DAB	267515	1
Bombastus	Zinc paste, soft DAB	280234	1
Bombastus	Zinc paste, soft DAB	278228	2
Bombastus	Zinc paste, soft DAB	270463	3
Bombastus	Zinc paste, soft DAB	276540	3
Bombastus	Zinc paste, soft DAB	271195	2
Bombastus	Zinc paste, soft DAB	271194	1
Bombastus	Zinc paste, soft DAB	267591	1
Bombastus	Zinc paste, soft DAB	266383	2
Bombastus	Zinc paste, soft DAB	11269801	1
Bombastus	Zinc paste, soft DAB	271343	2
Bombastus	Zinc paste, soft DAB	274971	1
Bombastus	Zinc paste, soft DAB	273701	1
Bombastus	Zinc paste, soft DAB	281117	1
Bombastus	Zinc paste, soft DAB	264706	1
Bombastus	Zinc paste, soft DAB	262918	1
Bombastus	Zinc paste, soft DAB	274972	2
Bombastus	Zinc paste, soft DAB	273704	3
Brückenapotheke	Hydrophilic metronida...	3913E-07001	1
Brückenapotheke	Hydrophilic metronida...	5013A-07001	1
Brückenapotheke	Hydrophilic metronida...	414M-07001	1
Brückenapotheke	Hydrophilic metronida...	3113Q-07001	1
Brückenapotheke	Hydrophilic metronida...	3313M-07001	1
Brücken-Apotheke	Hydrophilic metronida...	4613M-07001	1
Caelo	Acorus tincture (70%)...	19121322	1
Caelo	Acorus tincture (70%)...	11355115	1
Caelo	Alcoholic liniment	12103312	1
Caelo	Alcoholic liniment	34365645645	1
Caelo	Alcoholic liniment	445255	1
Caelo	Alcoholic liniment	5213E2316	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Alcoholic liniment	12103319	2
Caelo	Alcoholic liniment	2612M-02925	1
Caelo	Alcoholic liniment	204312316	1
Caelo	Alcoholic liniment	2507212316	1
Caelo	Aloe vera oil	11014807	2
Caelo	Aloe vera oil	11014806	1
Caelo	Aloe vera oil	11014804	1
Caelo	Aloe vera oil	12033107	1
Caelo	Aloe vera oil	12033114	2
Caelo	Aloe vera oil	12033111	1
Caelo	Aloe vera oil	123033108	1
Caelo	Aloe vera oil	12383702	2
Caelo	Aloe vera oil	12383703	2
Caelo	Aloe vera oil	12033101	2
Caelo	Aloe vera oil	12033104	1
Caelo	Aloe vera oil	12383706	5
Caelo	Aloe vera oil	12033113	4
Caelo	Aloe vera oil	12033106	1
Caelo	Aloe vera oil	12033108	1
Caelo	Apricot aroma	14C19-B30-292291	1
Caelo	Arnica oil, fat	11096702	1
Caelo	Arnica oil, fat	Sanacorb	1
Caelo	Arnica oil, fat	13186602	2
Caelo	Arnica oil, fat	12088504	3
Caelo	Arnica oil, fat	12359703	1
Caelo	Arnica oil, fat	13042503	3
Caelo	Arnica oil, fat	13186601	2
Caelo	Arnica oil, fat	13186604	1
Caelo	Arnica tincture (70%)...	12301009	1
Caelo	Arnica tincture (70%)...	13094613	1
Caelo	Arnica tincture (70%)...	10362809	1
Caelo	Arnica tincture (70%)...	20090312-4	1
Caelo	Arnica tincture (70%)...	1213A-02398	1
Caelo	Arnica tincture (70%)...	12031017	8
Caelo	Arnica tincture (70%)...	13094610	2
Caelo	Arnica tincture (70%)...	13094609	5
Caelo	Arnica tincture (70%)...	26031301	1
Caelo	Arnica tincture (70%)...	13094607	3
Caelo	Arnica tincture (70%)...	10341801	1
Caelo	Arnica tincture (70%)...	13094612	2
Caelo	Arnica tincture (70%)...	12031004	2
Caelo	Arnica tincture (70%)...	10362815	2
Caelo	Arnica tincture (70%)...	12031009	5
Caelo	Arnica tincture (70%)...	10362814	3
Caelo	Arnica tincture (70%)...	10362805	1
Caelo	Arnica tincture (70%)...	8111301	1
Caelo	Arnica tincture (70%)...	10362820	1
Caelo	Arnica tincture (70%)...	12031001	1
Caelo	Arnica tincture (70%)...	12031008	1
Caelo	Arnica tincture (70%)...	12031015	7
Caelo	Basil oil	13153601	1
Caelo	Belladonna tincture, ...	12315111	2
Caelo	Belladonna tincture, ...	12315104	2
Caelo	Belladonna tincture, ...	11322601	1
Caelo	Belladonna tincture, ...	11322611	2
Caelo	Belladonna tincture, ...	13273004	1
Caelo	Belladonna tincture, ...	11322615	1
Caelo	Belladonna tincture, ...	12315101	1
Caelo	Benzalkonium chloride	13307902	1
Caelo	Benzalkonium chloride	12033411	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Benzalkonium chloride	13070201	1
Caelo	Benzalkonium chloride	13070206	6
Caelo	Benzalkonium chloride	12033403	2
Caelo	Benzoic lard (stabili...	12063902	1
Caelo	Bitter orange tinctur...	12231906	2
Caelo	Bitter orange tinctur...	27072013	1
Caelo	Bitter orange tinctur...	13003301	1
Caelo	Bitter orange tinctur...	13003309	1
Caelo	Bitter orange tinctur...	13003302	1
Caelo	Calendula mother tinc...	13142811	6
Caelo	Calendula mother tinc...	12117313	1
Caelo	Calendula oil, fat	13185401	2
Caelo	Calendula oil, fat	nachprüfung	1
Caelo	Calendula oil, fat	13069603	1
Caelo	Calendula oil, fat	12373101	2
Caelo	Calendula oil, fat	12373103	1
Caelo	Calendula oil, fat	12308502	2
Caelo	Calendula oil, fat	13251002	3
Caelo	Calendula oil, fat	11183102	1
Caelo	Calendula oil, fat	12216904	1
Caelo	Calendula oil, fat	12189901	1
Caelo	Calendula oil, fat	12373102	3
Caelo	Calendula oil, fat	12216901	2
Caelo	Calendula oil, fat	12189903	3
Caelo	Calendula oil, fat	13185402	4
Caelo	Calendula oil, fat	13251002/1	1
Caelo	Calendula oil, fat	13251003	4
Caelo	Calendula oil, fat	13355701	1
Caelo	Calendula oil, fat	11274801	2
Caelo	Calendula oil, fat	12006802	4
Caelo	Calendula oil, fat	13069601	1
Caelo	Calendula oil, fat	14001004	2
Caelo	Calendula oil, fat	122116904	1
Caelo	Calendula oil, fat	13355702	2
Caelo	Calendula tincture (7...	11332702	1
Caelo	Calendula tincture (7...	11143803	1
Caelo	Calendula tincture (7...	12117301	3
Caelo	Calendula tincture (7...	13142804	7
Caelo	Calendula tincture (7...	12117305	1
Caelo	Calendula tincture (7...	12117308	1
Caelo	Calendula tincture (7...	12117310	2
Caelo	Calendula tincture (7...	12117311	7
Caelo	Calendula tincture (7...	12117320	1
Caelo	Calendula tincture (7...	12117317	5
Caelo	Calendula tincture (7...	12117312	3
Caelo	Calendula tincture (7...	12117313	1
Caelo	Calendula tincture (7...	13142808	12
Caelo	Calendula tincture (7...	19061311	1
Caelo	Calendula tincture (7...	11332710	1
Caelo	Calendula tincture (7...	10260111	1
Caelo	Calendula tincture (7...	13142802	2
Caelo	Camphor oil 10%	15028004	2
Caelo	Camphor oil 10%	15338202	2
Caelo	Camphor oil 10%	12269501	1
Caelo	Camphor oil 10%	13317401	4
Caelo	Camphor oil 10%	13043901	3
Caelo	Camphor oil 10%	18082015A	1
Caelo	Camphor oil 10%	15338202	1
Caelo	Camphor oil 10%	15338201	1
Caelo	Camphor oil 10%	15208001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Camphor oil 10%	15028002	1
Caelo	Camphor oil 10%	14139201	1
Caelo	Camphor oil 10%	15028003	1
Caelo	Camphor oil 10%	14139202	2
Caelo	Camphor oil 10%	14139202	3
Caelo	Camphor oil 10%	13317402	10
Caelo	Camphor oil 10%	17003105	1
Caelo	Camphor oil 10%	16085202	1
Caelo	Camphor oil 10%	15338205	1
Caelo	Camphor oil 10%	15028003	1
Caelo	Carbomer gel 1%	1102D002	1
Caelo	Carotene (ca. 4m IU)	13008404	1
Caelo	Castor oil, hydrated	12114602	3
Caelo	Castor oil, hydrated	16071201	1
Caelo	Castor oil, hydrated	11359502	2
Caelo	Castor oil, hydrated	14199604	1
Caelo	Castor oil, hydrated	141200-02	1
Caelo	Castor oil, hydrated	1304013-01	1
Caelo	Castor oil, hydrated	12114606	1
Caelo	Castor oil, hydrated	12114605	2
Caelo	Castor oil, hydrated	12071101	1
Caelo	Castor oil, hydrated	12042103	2
Caelo	Castor oil, hydrated	11359508	6
Caelo	Castor oil, hydrated	11359507	1
Caelo	Castor oil, hydrated	11260603	2
Caelo	Castor oil, hydrated	11260609	3
Caelo	Castor oil, hydrated	11260604	1
Caelo	Castor oil, hydrated	1112002-01	1
Caelo	Castor oil, hydrated	8061201	1
Caelo	Castor oil, hydrated	4101201	1
Caelo	Castor oil, hydrated	291722	1
Caelo	Castor oil, hydrated	11359501	2
Caelo	Castor oil, hydrated	25011301	1
Caelo	Cayenne pepper tinctu...	12207111	2
Caelo	Cayenne pepper tinctu...	11102425	1
Caelo	Cayenne pepper tinctu...	11102419	1
Caelo	Cayenne pepper tinctu...	13072013	1
Caelo	Cinchona bark tinctur...	13069808	1
Caelo	Coal tar solution	15203703	1
Caelo	Coal tar solution	15203704	1
Caelo	Coal tar solution	24131-01799	1
Caelo	Coal tar solution	90391069	2
Caelo	Coal tar solution	nachprüfung	1
Caelo	Coal tar solution	9345406	1
Caelo	Coal tar solution	9345509	2
Caelo	Coal tar solution	9345510	1
Caelo	Coal tar solution	27071207	1
Caelo	Coal tar solution	14331509	1
Caelo	Coal tar solution	14331403	1
Caelo	Coal tar solution	14331402	1
Caelo	Coal tar solution	14331401	1
Caelo	Coal tar solution	13193713	2
Caelo	Coal tar solution	13193703	4
Caelo	Coal tar solution	11319209	1
Caelo	Coal tar solution	11319207	12
Caelo	Coal tar solution	11319203	1
Caelo	Coal tar solution	113192017	1
Caelo	Coal tar solution	9345507	1
Caelo	Coal tar solution	9039201	1
Caelo	Coal tar solution	12286103	20

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Coal tar solution	14331506	1
Caelo	Coconut oil, hardened	11379517	1
Caelo	Coconut oil, hardened	11379511	1
Caelo	Coconut oil, hardened	11379510	2
Caelo	Coconut oil, hardened	13010804	1
Caelo	Coconut oil, hardened	13010808	2
Caelo	Coconut oil, hardened	613E-07004	1
Caelo	Comfrey root tincture...	12054012	2
Caelo	Comfrey root tincture...	12054008	1
Caelo	Comfrey root tincture...	11121314	2
Caelo	Comfrey tincture, cle...	4512E-02482	1
Caelo	Comfrey tincture, cle...	12177005	1
Caelo	Comfrey tincture, cle...	14298513	1
Caelo	Comfrey tincture, cle...	11121406	1
Caelo	Comfrey tincture, cle...	12177001	3
Caelo	Food colorant blue (...)	12322306	1
Caelo	Food colorant blue (...)	12038601	1
Caelo	Food colorant blue (...)	13212405	1
Caelo	Food colorant green ...	13020105	1
Caelo	Food colorant red (E...	13020001	1
Caelo	Food colorant red (E...	13020002	1
Caelo	Food colorant red (E...	13020005	1
Caelo	Food colorant yellow...	12022212	1
Caelo	Food colorant yellow...	13075706	1
Caelo	Formic acid 85%	13373301	1
Caelo	Formic acid 85%	10343801	1
Caelo	Gentian tincture (70%...	12035507	4
Caelo	Gentian tincture (70%...	10334609	2
Caelo	Gentian tincture (70%...	11220504	1
Caelo	Gentian tincture (70%...	13036405	4
Caelo	Gentian tincture (70%...	300711	1
Caelo	Gentian tincture (70%...	4109M-02442	1
Caelo	Gentian tincture (70%...	11220506	1
Caelo	Gentian tincture (70%...	13036409	2
Caelo	Ginger tincture (70%)...	11323810	1
Caelo	Ginger tincture (70%)...	11053812	1
Caelo	Hawthorn tincture of ...	12211008	1
Caelo	Horse chestnut seed l...	12293604	1
Caelo	Horse chestnut seed l...	11194611	2
Caelo	Horse chestnut seed t...	12295308	1
Caelo	Hydrogen peroxide sol...	7575	1
Caelo	Hydrogen peroxide sol...	12091302	1
Caelo	Hydrogen peroxide sol...	1003	1
Caelo	Hydrogen peroxide sol...	610093	1
Caelo	Hydrogen peroxide sol...		1
Caelo	Hydrophobic triclosan...	S011503.1	1
Caelo	Hydrophobic triclosan...	S0211201.1	1
Caelo	Hydrophobic triclosan...	3131-07001	1
Caelo	Hydrophobic triclosan...	2111310005	1
Caelo	Hydroxyethylcellulose...	14081202	1
Caelo	Hydroxyethylcellulose...	211111304	1
Caelo	Hydroxyethylcellulose...	8101304	1
Caelo	Hydroxyethylcellulose...	13084404	13
Caelo	Hydroxyethylcellulose...	11341501	1
Caelo	Hydroxyethylcellulose...	8101303	1
Caelo	Hydroxyethylcellulose...	6121303	1
Caelo	Hydroxyethylcellulose...	6061308	1
Caelo	Hydroxyethylcellulose...	13009303	3
Caelo	Hydroxyethylcellulose...	12351501	5
Caelo	Hydroxyethylcellulose...	13084403	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydroxyethylcellulose...	13084402	12
Caelo	Hydroxyethylcellulose...	13021302	1
Caelo	Hydroxyethylcellulose...	12351502	4
Caelo	Hydroxyethylcellulose...	13009501	10
Caelo	Hydroxyethylcellulose...	6121301	1
Caelo	Hydroxyethylcellulose...	13009301	5
Caelo	Hydroxyethylcellulose...	5121205	1
Caelo	Hydroxyethylcellulose...	6051303	1
Caelo	Hydroxyethylcellulose...	14066401	2
Caelo	Hydroxyethylcellulose...	14066301	3
Caelo	Hydroxyethylcellulose...	13393301	1
Caelo	Hydroxyethylcellulose...	14111307	1
Caelo	Hydroxyethylcellulose...	14111308	1
Caelo	Hydroxyethylcellulose...	13393201	1
Caelo	Hydroxyethylcellulose...	13286801	2
Caelo	Hydroxyethylcellulose...	17101207	1
Caelo	Hydroxyethylcellulose...	13286601	5
Caelo	Hydroxyethylcellulose...	13286401	1
Caelo	Hydroxyethylcellulose...	13286301	3
Caelo	Hydroxyethylcellulose...	6061307	1
Caelo	Hydroxyethylcellulose...	13286201	5
Caelo	Hydroxyethylcellulose...	1319401	1
Caelo	Hydroxyethylcellulose...	13190901	3
Caelo	Hydroxyethylcellulose...	13190801	5
Caelo	Hydroxyethylcellulose...	13190701	4
Caelo	Hydroxyethylcellulose...	17101201	1
Caelo	Hydroxyethylcellulose...	2101201	1
Caelo	Hydroxyethylcellulose...	2101203	1
Caelo	Hydroxyethylcellulose...	3286101	1
Caelo	Hydroxyethylcellulose...	3786917	1
Caelo	Hydroxyethylcellulose...	4091303	1
Caelo	Hydroxyethylcellulose...	10111203	1
Caelo	Hydroxyethylcellulose...	13286101	6
Caelo	Hydroxyethylcellulose...	21111305	1
Caelo	Hydroxyethylcellulose...	10111202	1
Caelo	Hydroxyethylcellulose...	11011402	1
Caelo	Hydroxyethylcellulose...	1219401	1
Caelo	Hydroxyethylcellulose...	12209201	1
Caelo	Hydroxyethylcellulose...	14200101	2
Caelo	Hydroxyethylcellulose...	12219201	6
Caelo	Hydroxyethylcellulose...	12219302	4
Caelo	Hydroxyethylcellulose...	12219501	5
Caelo	Hydroxyethylcellulose...	12244002	5
Caelo	Hydroxyethylcellulose...	30071302	1
Caelo	Hydroxyethylcellulose...	30071301	1
Caelo	Hydroxyethylcellulose...	29061302	1
Caelo	Hydroxyethylcellulose...	29061301	1
Caelo	Hydroxyethylcellulose...	29011301	1
Caelo	Hydroxyethylcellulose...	12320601	5
Caelo	Hydroxyethylcellulose...	12128202	2
Caelo	Hydroxyethylcellulose...	15067801	1
Caelo	Hydroxyethylcellulose...	26031306	1
Caelo	Hydroxyethylcellulose...	25091209	1
Caelo	Hydroxyethylcellulose...	20071204	1
Caelo	Hydroxyethylcellulose...	24011404	1
Caelo	Hydroxyethylcellulose...	24011403	1
Caelo	Hydroxyethylcellulose...	24011402	1
Caelo	Hydroxyethylcellulose...	23041305	1
Caelo	Hydroxyethylcellulose...	15067901	2
Caelo	Hydroxyethylcellulose...	16127002	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydroxyethylcellulose...	22061204	1
Caelo	Hydroxyethylcellulose...	15101311	1
Caelo	Hydroxyethylcellulose...	16061204	1
Caelo	Hydroxyethylcellulose...	21121201	1
Caelo	Hydroxyethylcellulose...	26101205	1
Caelo	Hydroxyethylcellulose...	12128201	5
Caelo	Hydroxyethylcellulose...	12128001	4
Caelo	Hydroxyethylcellulose...	12071701	2
Caelo	Hydroxyethylcellulose...	11264501	1
Caelo	Hydroxyethylcellulose...	11341701	1
Caelo	Hydroxyethylcellulose...	11341801	4
Caelo	Hydroxyethylcellulose...	11363302	1
Caelo	Hydroxyethylcellulose...	17101206	1
Caelo	Hydroxyethylcellulose...	11363401	1
Caelo	Hydroxyethylcellulose...	14081203	1
Caelo	Hydroxyethylcellulose...	12071501	3
Caelo	Hydroxyethylcellulose...	12071401	2
Caelo	Iodine tincture	28061302	1
Caelo	Iodine tincture	24011405	1
Caelo	Iodine tincture	14111310	1
Caelo	Iodine tincture	14111309	1
Caelo	Iodine tincture	13193503	1
Caelo	Iodine tincture	13193502	1
Caelo	Iodine tincture	9356601	1
Caelo	Iodine tincture	12285906	3
Caelo	Iodine tincture	13193501	2
Caelo	Iodine tincture	12285902	1
Caelo	Iodine tincture	11303104	1
Caelo	Iodine tincture	5021402	1
Caelo	Iodine tincture	5021308	1
Caelo	Iodine tincture	2101302	1
Caelo	Iodine tincture	1071303	1
Caelo	Iodine tincture	12285901	1
Caelo	Iodine tincture	12285903	1
Caelo	Lard	11263706	1
Caelo	Lard	11263707	1
Caelo	Lard	11263715	1
Caelo	Lard	11263714	1
Caelo	Lard	11379605	2
Caelo	Lard	11379606	1
Caelo	Lavender spirit	13166502	1
Caelo	Lavender spirit	13300803	2
Caelo	Lime water	13028602	1
Caelo	Macadamia oil, refined	14179304	2
Caelo	Macadamia oil, refined	13012106	1
Caelo	Macadamia oil, refined	13012104	1
Caelo	Marian thistle seed t...	13044304	1
Caelo	Marian thistle seed t...	12043822	1
Caelo	Marian thistle seed t...	11061110	1
Caelo	Matricaria flower tin...	3051310	1
Caelo	Matricaria flower tin...	3081201	1
Caelo	Matricaria flower tin...	9061201	1
Caelo	Matricaria flower tin...	11011303	1
Caelo	Matricaria flower tin...	11011304	1
Caelo	Matricaria flower tin...	11011307	1
Caelo	Matricaria flower tin...	13044206	1
Caelo	Matricaria flower tin...	11269905	1
Caelo	Matricaria flower tin...	11269911	2
Caelo	Matricaria flower tin...	24091208	1
Caelo	Matricaria flower tin...	1794290	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Matricaria flower tin...	11269917	3
Caelo	Matricaria flower tin...	11269920	5
Caelo	Matricaria flower tin...	11269924	1
Caelo	Matricaria flower tin...	14078402	1
Caelo	Matricaria flower tin...	13142803	1
Caelo	Matricaria flower tin...	12111305	1
Caelo	Matricaria flower tin...	13044202	1
Caelo	Matricaria liquid ext...	12360310	3
Caelo	Matricaria liquid ext...	10382713	1
Caelo	Matricaria liquid ext...	10305319	1
Caelo	Matricaria liquid ext...	1382716	1
Caelo	Matricaria liquid ext...	12360306	1
Caelo	Matricaria liquid ext...	12360303	2
Caelo	Matricaria liquid ext...	11303606	2
Caelo	Matricaria liquid ext...	10382716	3
Caelo	Matricaria liquid ext...	12254801	2
Caelo	Matricaria liquid ext...	12254803	1
Caelo	Matricaria oil, fat	12372903	5
Caelo	Matricaria oil, fat	11261506	1
Caelo	Matricaria oil, fat	1218903	1
Caelo	Melissa leaf tincture...	1206214	1
Caelo	Melissa leaf tincture...	3051307	1
Caelo	Melissa leaf tincture...	5121303	1
Caelo	Melissa leaf tincture...	11011305	1
Caelo	Melissa leaf tincture...	12206201	3
Caelo	Melissa leaf tincture...	10105412	1
Caelo	Melissa leaf tincture...	12206205	4
Caelo	Melissa leaf tincture...	11102812	1
Caelo	Melissa leaf tincture...	12206218	5
Caelo	Melissa leaf tincture...	12233609	2
Caelo	Melissa leaf tincture...	12233612	1
Caelo	Melissa leaf tincture...	1223612	1
Caelo	Melissa leaf tincture...	1233609	1
Caelo	Melissa leaf tincture...	27111307	1
Caelo	Melissa leaf tincture...	12206214	6
Caelo	Myrrh tincture (90%) ...	20111214	1
Caelo	Myrrh tincture (90%) ...	1910211095	1
Caelo	Myrrh tincture (90%) ...	24081209	1
Caelo	Myrrh tincture (90%) ...	3101212460	1
Caelo	Myrrh tincture (90%) ...	10112306	1
Caelo	Myrrh tincture (90%) ...	614M2460	1
Caelo	Myrrh tincture (90%) ...	13273301	7
Caelo	Myrrh tincture (90%) ...	2704212460	1
Caelo	Myrrh tincture (90%) ...	11260918	3
Caelo	Myrrh tincture (90%) ...	11260908	1
Caelo	Myrrh tincture (90%) ...	12197106	1
Caelo	Myrrh tincture (90%) ...	45132U2460	2
Caelo	Myrrh tincture (90%) ...	2913M2460	2
Caelo	Myrrh tincture (90%) ...	230312	1
Caelo	Myrrh tincture (90%) ...	20101012-1	1
Caelo	Myrrh tincture (90%) ...	2002312460	2
Caelo	Myrrh tincture (90%) ...	15110504	1
Caelo	Myrrh tincture (90%) ...	14244206	1
Caelo	Myrrh tincture (90%) ...	133273301	1
Caelo	Myrrh tincture (90%) ...	13273308	4
Caelo	Myrrh tincture (90%) ...	13273303	4
Caelo	Myrrh tincture (90%) ...	13273302	9
Caelo	Myrrh tincture (90%) ...	12197305	8
Caelo	Myrrh tincture (90%) ...	12197303	16
Caelo	Myrrh tincture (90%) ...	12197204	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Myrrh tincture (90%) ...	12197107	2
Caelo	Myrrh tincture (90%) ...	12197104	2
Caelo	Myrrh tincture (90%) ...	12197103	2
Caelo	Myrrh tincture (90%) ...	11260909	3
Caelo	Myrrh tincture (90%) ...	11260914	5
Caelo	Myrrh tincture (90%) ...	12197105	12
Caelo	Myrrh tincture (90%) ...	11260916	2
Caelo	Myrrh tincture (90%) ...	11260910	1
Caelo	Myrrh tincture (90%) ...	12197112	3
Caelo	Myrrh tincture (90%) ...	11019810	1
Caelo	Myrrh tincture (90%) ...	13273404	1
Caelo	Myrrh tincture (90%) ...	nachprüfung	1
Caelo	Myrrh tincture (90%) ...	11019809	1
Caelo	Nettle root tincture ...	11187409	1
Caelo	Nettle root tincture ...	12341302	1
Caelo	Nettle tincture (70%)...	11108810	1
Caelo	Nettle tincture (70%)...	1108810	1
Caelo	Opium tincture, stand. ...	13144005	3
Caelo	Opium tincture, stand. ...	12380307	1
Caelo	Opium tincture, stand. ...	12380302	1
Caelo	Opium tincture, stand. ...	12162608	1
Caelo	Opium tincture, stand. ...	12162602	1
Caelo	Opium tincture, stand. ...	12121305	1
Caelo	Opium tincture, stand. ...	10011404	1
Caelo	Opium tincture, stand. ...	8061301	1
Caelo	Opium tincture, stand. ...	13144010	1
Caelo	Opium tincture, stand. ...	6009916	1
Caelo	Opium tincture, stand. ...	5008715	3
Caelo	Opium tincture, stand. ...	5005712	8
Caelo	Opium tincture, stand. ...	4007614	4
Caelo	Opium tincture, stand. ...	3071303	1
Caelo	Opium tincture, stand. ...	2007414	1
Caelo	Opium tincture, stand. ...	2007314	1
Caelo	Opium tincture, stand. ...	8061302	1
Caelo	Opium tincture, stand. ...	13144011	1
Caelo	Opium tincture, stand. ...	13144004	1
Caelo	Opium tincture, stand. ...	13394902	1
Caelo	Opium tincture, stand. ...	11007213	10
Caelo	Opium tincture, stand. ...	30031301	1
Caelo	Opium tincture, stand. ...	13260602	1
Caelo	Opium tincture, stand. ...	25071302	1
Caelo	Opium tincture, stand. ...	25071303	1
Caelo	Opium tincture, stand. ...	14108510	1
Caelo	Opium tincture, stand. ...	14108504	1
Caelo	Opium tincture, stand. ...	16121304	1
Caelo	Opium tincture, stand. ...	22031302	1
Caelo	Opium tincture, stand. ...	24011406	1
Caelo	Orange flower oil, ar. ...	1123005	1
Caelo	Orange flower oil, ar. ...	nachprüfung	1
Caelo	Orange flower oil, ar. ...		1
Caelo	Orange flower oil, ar. ...	13121302	1
Caelo	Orange flower oil, ar. ...	13154804	1
Caelo	Orange flower oil, ar. ...	11232005	3
Caelo	Orange flower oil, ar. ...	13154801	1
Caelo	Orange flower water	13339301	1
Caelo	Orange flower water	14156804	1
Caelo	Orange flower water	12328602	1
Caelo	Orange flower water	12032102	1
Caelo	Orange flower water	11327301	2
Caelo	Orange flower water	12359606	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Orange flower water	10371902	1
Caelo	Orange flower water	1339301	1
Caelo	Orange flower water	13096605	3
Caelo	Orange flower water	12328604	3
Caelo	Orange flower water	1235606	1
Caelo	Orange flower water	13057302	2
Caelo	Orange flower water	13339302	2
Caelo	Peanut oil, hydrated	11266904	1
Caelo	Peanut oil, hydrated	13172801	2
Caelo	Peanut oil, hydrated	11300302	1
Caelo	Peanut oil, hydrated	12298204	2
Caelo	Peanut oil, hydrated	12298205	2
Caelo	Peanut oil, hydrated	12298206	4
Caelo	Peanut oil, hydrated	10206801	2
Caelo	Peanut oil, hydrated	12298203	1
Caelo	Peanut oil, hydrated	13172803	4
Caelo	Peppermint spirit 80.2%	12327202	2
Caelo	Peppermint spirit 80.2%	11185901	1
Caelo	Peppermint spirit 80.2%	11185902	1
Caelo	Peppermint spirit 80.2%	11339501	1
Caelo	Peppermint tincture (...)	24091206	1
Caelo	Peppermint tincture (...)	11207414	3
Caelo	Peppermint tincture (...)	12293409	1
Caelo	Peppermint tincture (...)	11207410	1
Caelo	Peppermint tincture (...)	10106507	1
Caelo	Peppermint tincture (...)	9089504	1
Caelo	Peppermint tincture (...)	9089503	1
Caelo	Peru balsam, original	13196609	1
Caelo	Peru balsam, original	131966609	1
Caelo	Peru balsam, original	12020602	3
Caelo	Peru balsam, original	11150808	3
Caelo	Peru balsam, original	13196609	1
Caelo	Peru balsam, original	12302411	5
Caelo	Peru balsam, original	12302413	1
Caelo	Peru balsam, original	13196606	1
Caelo	Peru balsam, original	13196606	1
Caelo	Peru balsam, original	11025902	1
Caelo	Peru balsam, original	12020611	4
Caelo	Peru balsam, original	Retest	1
Caelo	Plantain herb liquid ...	11332621	1
Caelo	Plantain herb liquid ...	13008703	2
Caelo	Plantain herb liquid ...	13008706	1
Caelo	Propolis liquid extra. ...	14377411	1
Caelo	Propolis liquid extra. ...	14377408	1
Caelo	Propolis liquid extra. ...	12086102	1
Caelo	Propolis liquid extra. ...	12086105	1
Caelo	Propolis liquid extra. ...	12086115	1
Caelo	Propolis liquid extra. ...	13347601	1
Caelo	Propolis liquid extra. ...	13347613	1
Caelo	Propolis liquid extra. ...	12086101	1
Caelo	Propolis liquid extra. ...	12086113	1
Caelo	Propolis tincture (96%)	12206322	1
Caelo	Propolis tincture (96%)	10305414	1
Caelo	Propolis tincture (96%)	11269521	1
Caelo	Propolis tincture (96%)	11269511	1
Caelo	Propolis tincture (96%)	11269516	1
Caelo	Ratanhia tincture (70. ...)	11204310	1
Caelo	Ratanhia tincture (70. ...)	CH:008/0913	1
Caelo	Ratanhia tincture (70. ...)	14101305	1
Caelo	Ratanhia tincture (70. ...)	13051206	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ratanhia tincture (70...	13051205	1
Caelo	Ratanhia tincture (70...	13051202	1
Caelo	Ratanhia tincture (70...	12255403	2
Caelo	Ratanhia tincture (70...	11292802	6
Caelo	Rose water	1021301	1
Caelo	Rose water	11355401	1
Caelo	Rose water	13306603	5
Caelo	Rose water	12123203	1
Caelo	Rose water	12288902	2
Caelo	Rose water	12185405	3
Caelo	Rose water	13137305	6
Caelo	Rose water	1234201	1
Caelo	Rose water	13343604	2
Caelo	Rose water	13137303	1
Caelo	Rose water	13137302	2
Caelo	Rose water	12342406	2
Caelo	Rose water	12342405	1
Caelo	Rose water	12342404	2
Caelo	Rose water	12342402	2
Caelo	Rose water	12288903	9
Caelo	Rose water	12299802	5
Caelo	Rose water	11286005	1
Caelo	Rose water	31101305	1
Caelo	Rose water	15253701	1
Caelo	Rose water	13306604	7
Caelo	Rose water	13255902	4
Caelo	Rose water	1317302	1
Caelo	Rose water	12342407	2
Caelo	Rose water	12342401	2
Caelo	Rose water	12299801	2
Caelo	Rosemary leaf tinctur...	9135005	1
Caelo	Rosemary leaf tinctur...	9135008	1
Caelo	Rosemary leaf tinctur...	140613	1
Caelo	Sage leaf liquid extr...	11372217	1
Caelo	Sage leaf liquid extr...	11372211	2
Caelo	Sage leaf liquid extr...	11372208	1
Caelo	Sage tincture (70%) 1...	12196607	2
Caelo	Sage tincture (70%) 1...	11220909	1
Caelo	Sage tincture (70%) 1...	10305312	1
Caelo	Sage tincture (70%) 1...	12196602	7
Caelo	Sage tincture (70%) 1...	14078111	1
Caelo	Sage tincture (70%) 1...	13125513	3
Caelo	Sage tincture (70%) 1...	13125509	1
Caelo	Sage tincture (70%) 1...	13125503	3
Caelo	Sage tincture (70%) 1...	1219661	1
Caelo	Sage tincture (70%) 1...	11220910	1
Caelo	Sage tincture (70%) 1...	1120910	1
Caelo	Sage tincture (70%) 1...	10305316	1
Caelo	Sage tincture (70%) 1...	13125505	1
Caelo	Sage tincture (70%) 1...	13125504	1
Caelo	Sage tincture (70%) 1...	12196611	7
Caelo	Sage tincture (70%) 1...	11220911	1
Caelo	Sage tincture (70%) 1...	1120912	1
Caelo	Sage tincture (70%) 1...	9345011	1
Caelo	Sage tincture (70%) 1...	15081804	1
Caelo	Sage tincture (70%) 1...	27111308	1
Caelo	Sage tincture (70%) 1...	4111M-02476	1
Caelo	Siam benzoin tincture...	1228119	1
Caelo	Siam benzoin tincture...	12228113	1
Caelo	Sodium lactate soluti...	131000701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Spearmint water	13184503	1
Caelo	Spearmint water	13184501	1
Caelo	St. John's wort oil	10380614	1
Caelo	St. John's wort oil	12076912	3
Caelo	St. John's wort oil	12076904	3
Caelo	St. John's wort oil	12076913	8
Caelo	St. John's wort oil	5442802	2
Caelo	St. John's wort oil	12293919	2
Caelo	St. John's wort oil	13065104	5
Caelo	St. John's wort oil	12293918	5
Caelo	St. John's wort oil	12293917	4
Caelo	St. John's wort oil	12292904	1
Caelo	St. John's wort oil	12293911	8
Caelo	St. John's wort oil	5346791	1
Caelo	St. John's wort oil	11372104	4
Caelo	St. John's wort oil	12193404	1
Caelo	St. John's wort oil	12293902	1
Caelo	St. John's wort oil	4959702	5
Caelo	St. John's wort oil	11205310	1
Caelo	St. John's wort oil	13065105	9
Caelo	St. John's wort oil	13065106	5
Caelo	St. John's wort oil	2071201	1
Caelo	St. John's wort oil	11372102	2
Caelo	St. John's wort oil	2061205	1
Caelo	St. John's wort oil	12293903	2
Caelo	St. John's wort oil	13155803	1
Caelo	St. John's wort oil	13252506	1
Caelo	St. John's wort oil	13389218	1
Caelo	St. John's wort oil	14130405	1
Caelo	St. John's wort oil	13155802	3
Caelo	St. John's wort oil	1501302	1
Caelo	St. John's wort oil	15358605	1
Caelo	St. John's wort oil	15358606	1
Caelo	St. John's wort oil	16037003	1
Caelo	St. John's wort oil	28091204	1
Caelo	St. John's wort oil	11080517	1
Caelo	St. John's wort oil	11205309	5
Caelo	St. John's wort oil	13155804	8
Caelo	St. John's wort tinct...	11195216	2
Caelo	St. John's wort tinct...	11195217	1
Caelo	St. John's wort tinct...	122377006	1
Caelo	St. John's wort tinct...	12237708	1
Caelo	St. John's wort tinct...	12237706	1
Caelo	St. John's wort tinct...	12237709	1
Caelo	St. John's wort tinct...	11091303	1
Caelo	St. John's wort tinct...	11195213	1
Caelo	Standardised cayenne ...	12053902	2
Caelo	Texapon [®] CS paste	16268909	2
Caelo	Texapon [®] CS paste	16268907	2
Caelo	Texapon [®] CS paste	16268904	1
Caelo	Texapon [®] CS paste	16268903	1
Caelo	Texapon [®] CS paste	1577804	1
Caelo	Texapon [®] CS paste	15177804	2
Caelo	Texapon [®] CS paste	12292007	1
Caelo	Texapon [®] CS paste	12292006	1
Caelo	Theriac	12067207	2
Caelo	Theriac	13187201	3
Caelo	Theriac	12067214	1
Caelo	Theriac	12067205	1
Caelo	Theriac	12067204	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Theriac	12067201	2
Caelo	Theriac	1121206	1
Caelo	Theriac	nachprüfung	1
Caelo	Theriac	11121203	1
Caelo	Thyme liquid extract	6838282	1
Caelo	Thyme liquid extract	1313M-01420	1
Caelo	Thyme liquid extract	100312	1
Caelo	Thyme liquid extract	11090508	1
Caelo	Thyme liquid extract	12072102	2
Caelo	Thyme liquid extract	13036502	2
Caelo	Thyme liquid extract	13036510	1
Caelo	Thyme liquid extract	13036516	1
Caelo	Thyme liquid extract	11090509	2
Caelo	Thyme tincture (70%) ...	18101304	1
Caelo	Thyme tincture (70%) ...	10379102	1
Caelo	Thyme tincture (70%) ...	10379108	1
Caelo	Thyme tincture (70%) ...	137918	1
Caelo	Thyme tincture (70%) ...	18101305	1
Caelo	Thyme tincture (70%) ...	10043113	1
Caelo	Thyme tincture (70%) ...	21031311	1
Caelo	Tormentil tincture (7...	13003003	1
Caelo	Tormentil tincture (7...	3713Q2484	2
Caelo	Tormentil tincture (7...	13263003	1
Caelo	Tormentil tincture (7...	11278307	2
Caelo	Tormentil tincture (7...	4212E-02484	1
Caelo	Tormentil tincture (7...	12255204	1
Caelo	Tormentil tincture (7...	11061207	1
Caelo	Tormentil tincture (7...	614M2484	1
Caelo	Tormentil tincture (7...	13003010	2
Caelo	Tormentil tincture (7...	11061210	2
Caelo	Tormentil tincture (7...	3101212484	1
Caelo	Tormentil tincture (7...	1910212484	2
Caelo	Tormentil tincture (7...	11278310	1
Caelo	Tormentil tincture (7...	12255201	1
Caelo	Tormentil tincture (7...	12255202	1
Caelo	Tormentil tincture (7...	13003002	2
Caelo	Tormentil tincture (7...	13003006	1
Caelo	Tormentil tincture (7...	13003008	1
Caelo	Tormentil tincture (7...	2103312484	2
Caelo	Tormentil tincture (7...	13003013	1
Caelo	Tormentil tincture (7...	13263007	1
Caelo	Tormentil tincture (7...	13003018	3
Caelo	Tormentil tincture (7...	13003012	6
Caelo	Tormentil tincture (7...	13003009	3
Caelo	Tormentil tincture (7...	11278309	2
Caelo	Tormentil tincture (7...	13263004	1
Caelo	Tretinoin solution 3...	13C28-N01	12
Caelo	Tretinoin solution 3...	11K09-N-03	1
Caelo	Triamcinolone acetoni...	12301201	2
Caelo	Triamcinolone acetoni...	11329001	1
Caelo	Triamcinolone acetoni...	12250701	1
Caelo	Triamcinolone acetoni...	13184701	7
Caelo	Triamcinolone acetoni...	13275801	4
Caelo	Triamcinolone acetoni...	1332501	1
Caelo	Triamcinolone acetoni...	11258801	1
Caelo	Valerian tincture (70...	12113004	3
Caelo	Valerian tincture (70...	12196203	1
Caelo	Valerian tincture (70...	12196210	2
Caelo	Valerian tincture (70...	12294007	4
Caelo	Valerian tincture (70...	13003205	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Valerian tincture (70...	13094503	2
Caelo	Valerian tincture (70...	13094506	2
Caelo	Valerian tincture (70...	13094507	5
Caelo	Valerian tincture (70...	1309505	1
Caelo	Valerian tincture (70...	8011423	1
Caelo	Valerian tincture (70...	6099.99.99	1
Caelo	Valerian tincture (70...	11011308	1
Caelo	Valerian tincture (70...	11379204	2
Caelo	Valerian tincture (70...	12030904	1
Caelo	Valerian tincture (70...	12113002	1
Caelo	Valerian tincture (70...	12030902	2
Caelo	Valerian tincture (70...	11379201	3
Caelo	Valerian tincture (70...	11254309	1
Caelo	Valerian tincture (70...	11187103	1
Caelo	Valerian tincture (70...	13349203	2
Caelo	Valerian tincture (70...	13094505-3380319003	1
Caelo	Valerian tincture (70...	13052601	1
Caelo	Valerian tincture (70...	12196204	6
Caelo	Valerian tincture (70...	18101302	1
Caelo	Valerian tincture (70...	12030908	1
Caelo	Valerian tincture (70...	16030101	1
Caelo	Valerian tincture (70...	14353011	1
Caelo	Valerian tincture (70...	11379203	2
Caelo	Valerian tincture (70...	13094505	7
Caelo	Valerian tincture (70...	13003206	6
Caelo	Valerian tincture (70...	12294008	2
Caelo	Valerian tincture (70...	13184905	4
Caelo	Valerian tincture (70...	12294002	3
Caelo	Vitamin A-palmitate (...)	120404301	2
Caelo	Vitamin A-palmitate (...)	11178101	2
Caelo	Vitamin A-palmitate (...)	11178109	1
Caelo	Vitamin A-palmitate (...)	130577002	1
Caelo	Vitamin A-palmitate (...)	130570012132910	1
Caelo	Vitamin A-palmitate (...)	13057002	9
Caelo	Vitamin A-palmitate (...)	12132910	6
Caelo	Vitamin A-palmitate (...)	12132906	4
Caelo	Vitamin A-palmitate (...)	12132902	2
Caelo	Vitamin A-palmitate (...)	12132904	4
Caelo	Vitamin A-palmitate (...)	13296208	1
Caelo	Vitamin A-palmitate (...)	13057002	1
Caelo	Vitamin A-palmitate (...)	3815926	1
Caelo	Vitamin A-palmitate (...)	13296206	1
Caelo	Vitamin A-palmitate (...)	12132904	7
Caelo	Vitamin A-palmitate (...)	12132910	5
Caelo	Vitamin A-palmitate (...)	13057003	1
Caelo	Vitamin A-palmitate (...)	1013Q-02565	1
Caelo	Vitamin A-palmitate (...)	10316605	1
Caelo	Vitamin A-palmitate (...)	1357007	2
Caelo	Vitamin A-palmitate (...)	1117801	1
Caelo	Vitamin A-palmitate (...)	12132906	4
Caelo	Vitamin A-palmitate (...)	12132902	4
Caelo	Witch hazel bark liqu...	11165602	1
Caelo	Witch hazel bark liqu...	11165605	1
Caelo	Witch hazel bark liqu...	13044102	3
Caelo	Witch hazel bark liqu...	13044106	1
Caelo	Witch hazel bark liqu...	13044104	2
Caelo	Witch hazel bark water	11238906	3
Caelo	Witch hazel bark water	12376906	2
Caelo	Witch hazel bark water	12376907	9
Caelo	Witch hazel bark water	12376918	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Witch hazel bark water	12376905	2
Caelo	Witch hazel bark water	11238904	1
Caelo	Witch hazel bark water	15430206	1
Caelo	Witch hazel bark water	14093309	2
Caelo	Witch hazel bark water	14093305	2
Caelo	Witch hazel bark water	14093302	1
Caelo	Witch hazel bark water	13324411	1
Caelo	Witch hazel bark water	13324403	7
Caelo	Witch hazel bark water	13324401	1
Caelo	Witch hazel bark water	12376912	1
Caelo	Witch hazel bark water	12376909	3
Caelo	Witch hazel bark water	12359605	1
Caelo	Witch hazel bark water	12196908	3
Caelo	Witch hazel bark water	12196906	4
Caelo	Witch hazel bark water	12196905	3
Caelo	Witch hazel bark water	1219608	1
Caelo	Witch hazel bark water	11238910	1
Caelo	Witch hazel bark water	11062814	1
Caelo	Witch hazel bark water	12376912Caelo	1
Caelo	Witch hazel bark water	12376901	3
Caelo	Witch hazel leaf liqu...	11060505	1
Caelo	Witch hazel leaf liqu...	27081302	1
Caelo	Witch hazel leaf liqu...	31031404	1
Caelo	Witch hazel leaf liqu...	27081303	1
Caelo	Witch hazel leaf liqu...	16071401	1
Caelo	Witch hazel leaf liqu...	15397406	1
Caelo	Witch hazel leaf liqu...	14112302	1
Caelo	Witch hazel leaf liqu...	12165010	1
Caelo	Wormwood tincture (70...	92780419	1
Caelo	Wormwood tincture (70...	11145005	1
Caelo	Wormwood tincture (70...	11145004	1
Caelo	Wormwood tincture (70...	11145001	1
Caelo	Wormwood tincture (70...	12237907	1
Caelo	Wormwood tincture (70...	1114508	1
Caelo	Wormwood tincture (70...	11145009	1
Caelo	Wormwood tincture (70...	12237905	4
Caelo	Wormwood tincture (70...	11145008	2
Caelo	Wormwood tincture (70...	12237901	1
Caelo	Wormwood tincture (70...	10260210	2
Caelo	Wormwood tincture (70...	12237911	1
Caelo	Wormwood tincture (70...	13187404	1
Caelo	Wormwood tincture (70...	12237909	4
Caelo	Wormwood tincture (70...	11145013	3
Caelo	Zinc ointment DAB	1234602	1
Caelo	Zinc ointment DAB	12101302	1
Caelo	Zinc ointment DAB	20270912-5	1
Caelo	Zinc ointment DAB	12343602	8
Caelo	Zinc ointment DAB	13072702	1
Caelo	Zinc ointment DAB	7031303	1
Caelo	Zinc ointment DAB	11367501	3
Caelo	Zinc ointment DAB	12042502	1
Caelo	Zinc ointment DAB	16101202	1
Caelo	Zinc ointment DAB	13257003	1
Caelo	Zinc ointment DAB	13257001	10
Caelo	Zinc ointment DAB	13393903	1
Caelo	Zinc ointment DAB	12098402	4
Caelo	Zinc ointment DAB	12042503	3
Caelo	Zinc ointment DAB	15003802	1
Caelo	Zinc oxide cream 25% ...	279250	2
Caelo	Zinc oxide cream 25% ...	2795250	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Zinc oxide lotion, aq...	12307902	1
Caelo	Zinc oxide lotion, aq...	14275002	1
Caelo	Zinc oxide lotion, aq...	15003201	3
Caelo	Zinc oxide lotion, aq...	16003302	1
Caelo	Zinc oxide lotion, aq...	11184601	1
Caelo	Zinc oxide lotion, aq...	15061201	1
Caelo	Zinc oxide lotion, aq...	12309501	33
Caelo	Zinc oxide lotion, aq...	15154301	1
Caelo	Zinc oxide lotion, aq...	15273502	2
Caelo	Zinc oxide lotion, aq...	13018401	18
Caelo	Zinc oxide lotion, aq...	16003202	1
Caelo	Zinc oxide lotion, aq...	12211501	8
Caelo	Zinc oxide lotion, aq...	19041309	1
Caelo	Zinc oxide lotion, aq...	13018803	7
Caelo	Zinc oxide lotion, aq...	20122013F	1
Caelo	Zinc oxide lotion, aq...	13024401	6
Caelo	Zinc oxide lotion, aq...	38112A-01815	1
Caelo	Zinc oxide lotion, aq...	13282901	19
Caelo	Zinc oxide lotion, aq...	3101217039	1
Caelo	Zinc oxide lotion, aq...	12309802	10
Caelo	Zinc oxide lotion, aq...	28121204	1
Caelo	Zinc oxide lotion, aq...	132801	1
Caelo	Zinc oxide lotion, aq...	12309601	7
Caelo	Zinc oxide lotion, aq...	13087601	32
Caelo	Zinc oxide lotion, aq...	2713M9101	1
Caelo	Zinc oxide lotion, aq...	2713M1815	1
Caelo	Zinc oxide lotion, aq...	24122013A	1
Caelo	Zinc oxide lotion, aq...	24101311	1
Caelo	Zinc oxide lotion, aq...	20310812-3	1
Caelo	Zinc oxide lotion, aq...	13087701	20
Caelo	Zinc oxide lotion, aq...	13282801	6
Caelo	Zinc oxide lotion, aq...	13283103	1
Caelo	Zinc oxide lotion, aq...	4121204	1
Caelo	Zinc oxide lotion, aq...	711217039	1
Caelo	Zinc oxide lotion, aq...	10082601	1
Caelo	Zinc oxide lotion, aq...	1013I-01815	1
Caelo	Zinc oxide lotion, aq...	11071303	1
Caelo	Zinc oxide lotion, aq...	11071304	1
Caelo	Zinc oxide lotion, aq...	11184501	1
Caelo	Zinc oxide lotion, aq...	11281301	1
Caelo	Zinc oxide lotion, aq...	112814101	1
Caelo	Zinc oxide lotion, aq...	11329902	1
Caelo	Zinc oxide lotion, aq...	12065203	1
Caelo	Zinc oxide lotion, aq...	12065301	9
Caelo	Zinc oxide lotion, aq...	12065401	3
Caelo	Zinc oxide lotion, aq...	12111306	1
Caelo	Zinc oxide lotion, aq...	12197103	1
Caelo	Zinc oxide lotion, aq...	6061301	1
Caelo	Zinc oxide lotion, aq...	12211401	10
Caelo	Zinc oxide lotion, aq...	12211701	17
Caelo	Zinc oxide lotion, aq...	12211802	5
Caelo	Zinc oxide lotion, aq...	1230951	1
Caelo	Zinc oxide lotion, aq...	12309701	5
Caelo	Zinc oxide lotion, aq...	12309702	6
Caelo	Zinc oxide lotion, aq...	13018801	8
Caelo	Zinc oxide lotion, aq...	1313901	1
Caelo	Zinc oxide lotion, aq...	13139301	19
Caelo	Zinc oxide lotion, aq...	13139302	5
Caelo	Zinc oxide lotion, aq...	13139601	36
Caelo	Zinc oxide lotion, aq...	1319601	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Zinc oxide lotion, aq. ...	13282701	1
Caelo	Zinc oxide lotion, aq. ...	13283101	9
Caelo	Zinc oxide lotion, aq. ...	13283002	6
Caelo	Zinc oxide lotion, aq. ...	12211601	2
Caelo	Zinc oxide lotion, aq. ...	14010501	3
Caelo	Zinc oxide lotion, aq. ...	13283001	3
Caelo	Zinc oxide lotion, aq. ...	14116701	1
Caelo	Zinc oxide lotion, aq. ...	14010602	1
Caelo	Zinc oxide lotion, aq. ...	1109779	1
Caelo	Zinc oxide lotion, aq. ...	14010502	1
Caelo	Zinc oxide oil DAC	12351702	19
Caelo	Zinc oxide oil DAC	13196401	38
Caelo	Zinc oxide oil DAC	12289802	3
Caelo	Zinc oxide oil DAC	13307202	35
Caelo	Zinc oxide oil DAC	14032701	2
Caelo	Zinc oxide oil DAC	15170902	1
Caelo	Zinc oxide oil DAC	15173002	1
Caelo	Zinc oxide oil DAC	16093002	1
Caelo	Zinc oxide oil DAC	2612M-02026	1
Caelo	Zinc oxide oil DAC	3112M-02026	1
Caelo	Zinc oxide oil DAC	13196501	2
Caelo	Zinc oxide oil DAC	12004701	2
Caelo	Zinc oxide oil DAC	12160002	30
Caelo	Zinc oxide oil DAC	12160101	1
Caelo	Zinc oxide oil DAC	12289801	1
Caelo	Zinc oxide oil DAC	12004601	2
Caelo	Zinc oxide oil DAC	12128602	1
Caelo	Zinc oxide paste SR90	Test2	1
Caelo	Zinc oxide paste SR90	271195	1
Caelo	Zinc oxide paste SR90	244154	1
Caelo	Zinc oxide paste, sof. ...	12366801	1
Caelo	Zinc oxide shaking mi. ...	1324401	1
Caelo	Zinc oxide shaking mi. ...	13080201	3
Caelo	Zinc oxide shaking mi. ...	13244001	1
Caelo	Zinc oxide shaking mi. ...	1308021	1
Caelo	Zinc oxide shaking mi. ...	12261701	1
Caelo	Zinc paste DAB	14293401	1
Caelo	Zinc paste DAB	13201801	25
Caelo	Zinc paste DAB	14092301	1
Caelo	Zinc paste DAB	20271012-1	1
Caelo	Zinc paste DAB	13043301	4
Caelo	Zinc paste DAB	14011402	1
Caelo	Zinc paste DAB	10345201	3
Caelo	Zinc paste DAB	12271601	10
Caelo	Zinc paste DAB	12206501	11
Caelo	Zinc paste DAB	20021310	1
Caelo	Zinc paste DAB	11290901	1
Caelo	Zinc paste DAB	11321701	4
Caelo	Zinc paste DAB	13201201	1
Caelo	Zinc paste DAB	13201701	2
Caelo	Zinc paste DAB	903312054	1
Caelo	Zinc paste DAB	12262301	2
Caelo	Zinc paste DAB	10245901	2
Caelo	Zinc paste DAB	13345202	1
Caelo	Zinc paste DAB	12066102	1
Caelo	Zinc paste DAB	11364102	3
Caelo	Zinc paste DAB	14293302	1
Caelo	Zinc paste DAB	1321903	1
Caelo	Zinc paste DAB	16068501	1
Caelo	Zinc paste DAB	12147801	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Zinc paste DAB	14201201	1
Caelo	Zinc paste DAB	11290801	2
Caelo	Zinc paste DAB	16013302	1
Caelo	Zinc paste DAB	12206401	6
Caelo	Zinc paste DAB	31051302	1
Caelo	Zinc paste DAB	15121-02054	1
Caelo	Zinc paste DAB	12271702	12
Caelo	Zinc paste DAB	11321702	1
Caelo	Zinc paste, soft DAB	13203401	13
Caelo	Zinc paste, soft DAB	13203301	17
Caelo	Zinc paste, soft DAB	13203501	12
Caelo	Zinc paste, soft DAB	13036901	1
Caelo	Zinc paste, soft DAB	13089101	6
Caelo	Zinc paste, soft DAB	13089001	15
Caelo	Zinc paste, soft DAB	13088901	13
Caelo	Zinc paste, soft DAB	16042101	1
Caelo	Zinc paste, soft DAB	14023401	1
Caelo	Zinc paste, soft DAB	13316301	25
Caelo	Zinc paste, soft DAB	15073001	2
Caelo	Zinc paste, soft DAB	14296301	1
Caelo	Zinc paste, soft DAB	15005202	1
Caelo	Zinc paste, soft DAB	13026801	7
Caelo	Zinc paste, soft DAB	13040201	13
Caelo	Zinc paste, soft DAB	13027001	8
Caelo	Zinc paste, soft DAB	14111313	1
Caelo	Zinc paste, soft DAB	14173001	2
Caelo	Zinc paste, soft DAB	13088701	30
Caelo	Zinc paste, soft DAB	14173201	1
Caelo	Zinc paste, soft DAB	1417601	1
Caelo	Zinc paste, soft DAB	13388001	3
Caelo	Zinc paste, soft DAB	14181702	1
Caelo	Zinc paste, soft DAB	13088801	12
Caelo	Zinc paste, soft DAB	14296201	1
Caelo	Zinc paste, soft DAB	13316701	1
Caelo	Zinc paste, soft DAB	12366701	12
Caelo	Zinc paste, soft DAB	14296801	1
Caelo	Zinc paste, soft DAB	14117801	1
Caelo	Zinc paste, soft DAB	13026901	21
Caelo	Zinc paste, soft DAB	213E-02055	1
Caelo	Zinc paste, soft DAB	13316401	6
Caelo	Zinc paste, soft DAB	12068801	8
Caelo	Zinc paste, soft DAB	16121901	1
Caelo	Zinc paste, soft DAB	16001201	3
Caelo	Zinc paste, soft DAB	15072901	1
Caelo	Zinc paste, soft DAB	15072701	3
Caelo	Zinc paste, soft DAB	14172901	1
Caelo	Zinc paste, soft DAB	14023801	1
Caelo	Zinc paste, soft DAB	11338601	3
Caelo	Zinc paste, soft DAB	14023701	3
Caelo	Zinc paste, soft DAB	13387801	10
Caelo	Zinc paste, soft DAB	13316601	6
Caelo	Zinc paste, soft DAB	13203701	2
Caelo	Zinc paste, soft DAB	13088602	1
Caelo	Zinc paste, soft DAB	13088/01	1
Caelo	Zinc paste, soft DAB	13026702	5
Caelo	Zinc paste, soft DAB	12366801	12
Caelo	Zinc paste, soft DAB	14023601	1
Caelo	Zinc paste, soft DAB	11061301	1
Caelo	Zinc paste, soft DAB	5517500	5
Caelo	Zinc paste, soft DAB	5611900	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Zinc paste, soft DAB	13316501	7
Caelo	Zinc paste, soft DAB	1601312055	1
Caelo	Zinc paste, soft DAB	12366703	1
Caelo	Zinc paste, soft DAB		2
Caelo	Zinc paste, soft DAB	12271301	5
Caelo	Zinc paste, soft DAB	63075606	1
Caelo	Zinc paste, soft DAB	91107159	1
Caelo	Zinc paste, soft DAB	12206201	1
Caelo	Zinc paste, soft DAB	12121601	10
Caelo	Zinc paste, soft DAB	12121501	4
Caelo	Zinc paste, soft DAB	12068701	4
Caelo	Zinc paste, soft DAB	11338801	4
Caelo	Zinc paste, soft DAB	8081301	1
Caelo	Zinc paste, soft DAB	711212055	1
Caelo	Zinc paste, soft DAB	70913WZP	1
Caelo	Zinc paste, soft DAB	17111203	1
Caelo	Zinc paste, soft DAB	12069001	5
Caelo	Zinc paste, soft DAB	13004101	1
Caelo	Zinc paste, soft DAB	2713E2055	1
Caelo	Zinc paste, soft DAB	12207901	8
Caelo	Zinc paste, soft DAB	208212055	1
Caelo	Zinc paste, soft DAB	12271501	21
Caelo	Zinc paste, soft DAB	12121302	2
Caelo	Zinc paste, soft DAB	11338001	2
Caelo	Zinc paste, soft DAB	11269701	2
Caelo	Zinc paste, soft DAB	1212601	1
Caelo	Zinc paste, soft DAB	1403212055	1
Caelo	Zinc paste, soft DAB	2101212055	1
Caelo	Zinc paste, soft DAB	414U-02055	1
Caelo	Zinc paste, soft DAB	3413E-02055	1
Caelo	Zinc paste, soft DAB	12068901	2
Caelo	Zinc paste, soft DAB	1212401	1
Caelo	Zinc paste, soft DAB	11338401	2
Caelo	Zinc paste, soft DAB	6021305	1
Caelo	Zinc paste, soft DAB	12207801	1
Caelo	Zinc paste, soft DAB	12121401	7
Caelo	Zinc paste, soft DAB	12208101	6
Caelo	Zinc paste, soft DAB	12208201	17
Caelo	Zinc paste, soft DAB	11099003	1
Caelo	Zinc paste, soft DAB	12271401	11
Caelo	Zinc paste, soft DAB	1268901	1
Caelo	Zinc paste, soft DAB	12308307	1
Caelo	Zinc paste, soft DAB	28031401	1
Caelo	Zinc paste, soft DAB	11202101	1
Caelo	Zinc paste, soft DRF ...	12160302	1
Caelo	Zinc paste, soft DRF ...	13094201	2
Caelo	Zinc paste, soft DRF ...	12378402	1
Caesar & Loretz	Plantain herb liquid ...	13008704	2
Caesar & Loretz GmbH	Coal tar solution	72816068	1
Caesar & Loretz GmbH	Hydrophobic triclosan...	S041303.1	3
Caesar & Loretz GmbH	Matricaria oil, fat	13204304	2
Caesar & Loretz GmbH	Matricaria oil, fat	12189803	2
Caesar & Loretz GmbH	Vitamin A-palmitate (...)	10316605	1
Caesar & Loretz GmbH	Zinc oxide lotion, aq. ...	13405901	4
Caesar & Loretz GmbH	Zinc paste DAB	281972	1
Caesar & Loretz GmbH ...	Zinc paste, soft DAB	11269601	2
Caesar & Loretz GmbH, ...	Benzalkonium chloride	12033407	2
Caesar & Loretz/ Noweda	Hydroxyethylcellulose...	5B150213	1
Caesar & Loretz/ Noweda	Plantain herb liquid ...	11332611	1
Caesar & Loretz/ Noweda	Plantain herb liquid ...	13008708	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz/ Noweda	Plantain herb liquid ...	11332616	3
Caesar & Loretz/ Noweda	Plantain herb liquid ...	11332613	2
Caesar & Loretz/ Noweda	Plantain herb liquid ...	11332620	1
Caesar § Lorenz GmbH	Aloe vera oil	12383707	1
Caesar u Loretz	Calendula oil, fat	13251004	1
Caesar&Loretz GmbH/An...	Peppermint spirit 80.2%	13252802	1
Caesar&lorenz,Phoenix	Lard	11095310	1
Caesar/Lorenz	Propolis tincture (96%)	12206319	1
Caesar/Noweda	Gentian tincture (70%...	13036408	1
Caesar/Noweda	Ratanhia tincture (70...	13051210	1
Calo	Zinc oxide lotion, aq...	14211601	1
Cheplapharm Arzneimit...	Valerian tincture (70...	70713	1
Defektur	Hydrogen peroxide sol...	6726	1
Defektur	Hydrogen peroxide sol...	6935	1
Defektur	Hydrogen peroxide sol...	6629	1
Defektur	Hydrophilic metronida...	3212M-03387	1
Defektur	Hydrophilic metronida...	1813E-07001	1
Defektur	Hydrophilic metronida...	4912Q-07001	1
Defektur	Hydrophilic metronida...	4012A-07001	1
Defektur	Hydrophilic metronida...	3712Q-07001	1
Defektur	Hydrophilic metronida...	3512L-07001	1
Defektur	Hydrophilic metronida...	3412M-07001	1
Defektur	Hydrophilic metronida...	2513q-07001	1
Defektur	Hydrophilic metronida...	2713E-07001	1
Defektur	Hydrophilic metronida...	4112E-07001	1
Defektur	Hydrophilic metronida...	4612E-03387	1
Defektur	Hydrophilic metronida...	813A-07001	1
Defektur	Hydrophilic metronida...	4812E-07001	1
Defektur	Hydrophilic metronida...	1613a-07001	1
Defektur	Hydrophilic metronida...	1813Q-07001	1
Defektur	Hydrophilic metronida...	1913A-07001	1
Defektur	Oxytetracycline ointm...	55/12	1
Defektur	Oxytetracycline ointm...	50/13	1
Defektur	Oxytetracycline ointm...	12/13	1
Defektur	Oxytetracycline ointm...	49/12	1
Defektur	Oxytetracycline ointm...	43/13	1
Defektur	Oxytetracycline ointm...	41/13	1
Defektur	Oxytetracycline ointm...	8/13	1
Defektur	Oxytetracycline ointm...	25/13	1
Defektur	Triamcinolone acetoni...	4713A	1
Defektur	Triamcinolone acetoni...	4213E	1
Defektur	Triamcinolone acetoni...	3813Q	1
Defektur	Triamcinolone acetoni...	4913Q	1
Defektur	Zinc paste, soft DAB	32/13	1
Defektur	Zinc paste, soft DAB	11/13	1
Defektur	Zinc paste, soft DAB	56/13	1
Defektur	Zinc paste, soft DAB	54/12	1
dhu	Lavender spirit	1041012	4
Eigenherstellung	Hydrophilic metronida...	912Q-07001	1
Eigenherstellung	Opium tincture, stand...	70613PML01	2
Erika apotheke	Gentian tincture (70%...	1122504	1
Erika apotheke	Tormentil tincture (7...	13003001	1
Erika apotheke	Wormwood tincture (70...	12237906	1
Erika Apotheke	Orange flower water	12359601	1
Erika Apotheke	Wormwood tincture (70...	313E-07062	1
Erika Apotheke	Wormwood tincture (70...	12237902	3
EuRho	Castor oil, hydrated	1212008-01	1
EuRho	Valerian tincture (70...	1301003-01	1
EuRho/Phönix	Castor oil, hydrated	1412006-01	1
Euro OTC	Arnica tincture (70%)...	1112029-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Arnica tincture (70%)...	1306038-01	1
Euro OTC	Camphor oil 10%	15338204	2
Euro OTC	Castor oil, hydrated	1202041-02	1
Euro OTC	Castor oil, hydrated	120204-01	1
Euro OTC	Castor oil, hydrated	2507213553	1
Euro OTC	Castor oil, hydrated	1512212009	1
Euro OTC	Castor oil, hydrated	1109001-01	2
Euro OTC	Castor oil, hydrated	1412006-02	1
Euro OTC	Castor oil, hydrated	1910212009	1
Euro OTC	Castor oil, hydrated	1206039-01	2
Euro OTC	Castor oil, hydrated	1206039-02	2
Euro OTC	Castor oil, hydrated	1112002-02	1
Euro OTC	Coal tar solution	1110031-02	6
Euro OTC	Coal tar solution	1307027-01	2
Euro OTC	Coal tar solution	1110031-01	6
Euro OTC	Coal tar solution	1110031	2
Euro OTC	Sodium lactate soluti...	1310007-01	3
Euro OTC	Sodium lactate soluti...	1212004-01	2
Euro OTC	Tretinoin solution 3...	13c28-n01	1
Euro OTC	Valerian tincture (70...	1304041-01	1
Euro OTC	Vitamin A-palmitate (...)	1305003-01-342631115	1
Euro OTC	Vitamin A-palmitate (...)	1310030-01	1
Euro OTC	Vitamin A-palmitate (...)	1204043-01	3
Euro OTC	Vitamin A-palmitate (...)	1305003-01	3
Euro OTC	Vitamin A-palmitate (...)	1204043-01	5
Euro OTC	Zinc paste, soft DAB	121014-01	1
Euro OTC	Zinc paste, soft DAB	1108059-01	4
Euro OTC	Zinc paste, soft DAB	13211014-01	1
Euro OTC	Zinc paste, soft DAB	1310019-01	2
Euro OTC	Zinc paste, soft DAB	1211014-01	9
Euro OTC	Zinc paste, soft DAB	121101401	1
Euro OTC	Zinc paste, soft DAB	1009043-01	3
Euro OTC	Zinc paste, soft DAB	214I-02055	1
EUro/AHD	Vitamin A-palmitate (...)	1211032-01	2
Fagron	Carbomer gel 1%	1309D001	1
Fagron	Carbomer gel 1%	1302D002	7
Fagron	Carbomer gel 1%	1309D002	1
Fagron	Carbomer gel 1%	1211D010	3
Fagron	Castor oil, hydrated	10L02-N03	1
Fagron	Coal tar solution	4513E1799	1
Fagron	Coal tar solution	16102013B	1
Fagron	Coal tar solution	11/05-N12	1
Fagron	Coal tar solution	11-05-N12	1
Fagron	Coal tar solution	12E08-N08	1
Fagron	Coal tar solution	12F19-N08	5
Fagron	Coal tar solution	13H06-N02	11
Fagron	Coal tar solution	13H06N02	1
Fagron	Coal tar solution	18061301	1
Fagron	Coal tar solution	2011.07.0288	2
Fagron	Coal tar solution	11F21-N01	3
Fagron	Coal tar solution	20230712-6	1
Fagron	Coal tar solution	12L05-N02	4
Fagron	Coal tar solution	12H28-N05	1
Fagron	Coal tar solution	13I26-N02	1
Fagron	Hydrogen peroxide sol...	13G09-M07	1
Fagron	Hydrogen peroxide sol...	13G12-M06	1
Fagron	Hydrogen peroxide sol...	13G09-M03	1
Fagron	Hydrogen peroxide sol...	13D05-M04	1
Fagron	Hydrogen peroxide sol...	13G12-M09	1
Fagron	Hydrogen peroxide sol...	13A02-M10	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hydrogen peroxide sol. ...	13K07-M04	1
Fagron	Hydrogen peroxide sol. ...	12I06-M05	1
Fagron	Hydrogen peroxide sol. ...	13A02-M13	1
Fagron	Hydroxyethylcellulose...	13009401	3
Fagron	Hydroxyethylcellulose...	12351401	1
Fagron	Iodine tincture	11309402	1
Fagron	Methylrosanilinium ch. ...	12D19-N01	1
Fagron	Opium tincture, stand. ...	5007714	3
Fagron	Opium tincture, stand. ...	9006012/Maros	1
Fagron	Peru balsam, original	3913E-01162	1
Fagron	Peru balsam, original	12K30-N02	2
Fagron	Peru balsam, original	2215I-01162	1
Fagron	Peru balsam, original	12H29-N010432/549	1
Fagron	Peru balsam, original	12A20-N13	2
Fagron	Peru balsam, original	9G03-N09	1
Fagron	Peru balsam, original	11G22-N06	2
Fagron	Rose water	737724	1
Fagron	Rose water	13B01-N01	3
Fagron	Rose water	12H-20-N09	1
Fagron	Rose water	12H20-N09	4
Fagron	St. John's wort oil	5346701	7
Fagron	St. John's wort oil	12L17-N15	1
Fagron	St. John's wort oil	13A30-N09	1
Fagron	St. John's wort oil	12A30-N09	1
Fagron	Tretinoin solution 3. ...	13I11-N04	7
Fagron	Tretinoin solution 3. ...	1108D001	1
Fagron	Tretinoin solution 3. ...	1308D00213/11-N04	1
Fagron	Tretinoin solution 3. ...	12C08-N14	4
Fagron	Tretinoin solution 3. ...	1308D002	1
Fagron	Tretinoin solution 3. ...	1204D004	1
Fagron	Tretinoin solution 3. ...	13C25-N01	1
Fagron	Tretinoin solution 3. ...	1501D004	1
Fagron	Tretinoin solution 3. ...	19061306	1
Fagron	Tretinoin solution 3. ...	13/11-N04	3
Fagron	Tretinoin solution 3. ...	12G05-N15	7
Fagron	Triamcinolone acetoni. ...	214Q	1
Fagron	Vitamin A-palmitate (...)	13C19-N12	5
Fagron	Wormwood tincture (70. ...)	13187402	1
Fagron	Zinc oxide lotion, aq. ...	4213Q-01815	1
Fagron	Zinc oxide lotion, aq. ...	13K19-T02-014523	1
Fagron	Zinc oxide lotion, aq. ...	10H19-T01-000747	1
Fagron	Zinc oxide oil SR 90	n. Anbruch3Monate	1
Fagron	Zinc paste DAB	12I06-T01	1
Fagron	Zinc paste DAB	9k06-t02	1
Fagron	Zinc paste, soft DAB	12J26-T03	6
Fagron	Zinc paste, soft DAB	10D15-T01	1
Fagron	Zinc paste, soft DAB	12/26-T03	1
Fagron	Zinc paste, soft DAB	13H13-T05-006158	1
Fagron	Zinc paste, soft DAB	4113A-02055	1
Fagron	Zinc paste, soft DAB	12271502	5
Fiebig	Coal tar solution	9345501	3
Fiebig	Hydrogen peroxide sol. ...	61007B	7
Fiebig	Zinc paste DAB	268489	1
Fiebig 23.10.12 ek:6, ...	Theriac	12067202	1
Fiebig ek:5.69EUR 10.08. ...	Witch hazel bark water	11238911	2
Finze	Zinc oxide lotion, aq. ...	12065501	5
Fischar	Hydrogen peroxide sol. ...	6303042A	1
Fischar	Hydrogen peroxide sol. ...	631109	1
Fischar	Hydrogen peroxide sol. ...	6311092	2
Fischar	Hydrogen peroxide sol. ...	6312082	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischar	Hydrogen peroxide sol...	6312061	1
Fischar	Hydrogen peroxide sol...	6312077	2
Fischar	Hydrogen peroxide sol...	6312081	3
Fischar	Hydrogen peroxide sol...	6312073	2
Fischar	Hydrogen peroxide sol...	6313025	2
Fischar	Hydrogen peroxide sol...	6315084	1
Fischar	Hydrogen peroxide sol...	6314018	1
Fischar	Hydrogen peroxide sol...	6313061	5
Fischar/AHD	Hydrogen peroxide sol...	6313062	1
Fischar/Anzag	Hydrogen peroxide sol...	6313064	2
Fischar/G	Hydrogen peroxide sol...	6313033	1
Fischar/G	Hydrogen peroxide sol...	6313093	2
Fischar/G	Hydrogen peroxide sol...	6312072	1
Fischar/Gehe	Hydrogen peroxide sol...	6313102	1
Fischar/GEHE/16.04.2013	Hydrogen peroxide sol...	6312114	3
Fischar/GUMO /02.08.2...	Hydrogen peroxide sol...	6313027	6
Fischar/Sanacorp	Hydrogen peroxide sol...	6312111	1
Fischbar	Hydrogen peroxide sol...	4101202	1
Gehe	Bitter orange tinctur...	11187211	1
Gehe	Castor oil, hydrated	272591	4
Gehe	Coal tar solution	12A20-N04	1
Gehe	Coal tar solution	701972-0001	2
Gehe	Hydrogen peroxide sol...	6312074	4
Gehe	Hydrogen peroxide sol...	6314021	2
Gehe	Hydrophilic metronida...	2512Q-07001	2
Gehe	Hydrophilic metronida...	4813E-07001	1
Gehe	Hydrophilic metronida...	5513E-03387	1
Gehe	Hydrophilic metronida...	1013Q-07001	1
Gehe	Hydrophilic metronida...	114M-07001	1
Gehe	Hydrophilic metronida...	2813M-07001	1
Gehe	Hydrophilic metronida...	3613E-07001	1
Gehe	Hydrophilic metronida...	1213E-07001	2
Gehe	Hydrophilic metronida...	5112Q-07001	1
Gehe	Hydrophobic triclosan...	S031103.1	1
Gehe	Hydrophobic triclosan...	S091107.1	2
Gehe	Hydrophobic triclosan...	S101209.1	17
Gehe	Hydrophobic triclosan...	S021201.1	5
Gehe	Hydroxyethylcellulose...	22b070912	1
Gehe	Hydroxyethylcellulose...	5B090312	4
Gehe	Methoxsalen 0.006% Co...	13D004A	1
Gehe	Methoxsalen 0.006% Co...	13D004B	1
Gehe	Metronidazole 10% (no...	8B170712	1
Gehe	Metronidazole 10% (no...	4B170413	1
Gehe	Opium tincture, stand...	6007814	2
Gehe	Searing paste SR	4B080313	1
Gehe	Searing paste SR	10B030712	1
Gehe	Vitamin A-palmitate (...)	1211032-01	4
Gehe	Vitamin A-palmitate (...)	1110018-01	2
Gehe	Vitamin A-palmitate (...)	11178103	2
Gehe	Vitamin A-palmitate (...)	11178101	3
Gehe	Zinc oxide lotion, aq...	12211402	8
Gehe	Zinc oxide oil DAC	279563	3
Gehe	Zinc oxide oil DAC	282085	1
Gehe	Zinc oxide oil DAC	275085	1
Gehe	Zinc oxide oil DAC	289826	1
Gehe	Zinc oxide oil DAC	275086	2
Gehe	Zinc oxide paste SR90	264198	1
Gehe	Zinc oxide paste SR90	268240	1
Gehe	Zinc oxide paste SR90	266451	1
Gehe	Zinc oxide paste SR90	270919	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Zinc oxide paste SR90	279386	1
Gehe	Zinc oxide paste SR90	276199	2
Gehe	Zinc oxide paste, sof...	36B231012	1
Gehe	Zinc oxide paste, sof...	23B160713	1
Gehe	Zinc oxide paste, sof...	260912	1
Gehe	Zinc oxide paste, sof...	30071502	1
Gehe	Zinc paste DAB	271138	3
Gehe	Zinc paste DAB	276144	3
Gehe	Zinc paste, soft DAB	130112	1
Gehe	Zinc paste, soft DAB	12A250412	1
Gehe	Zinc paste, soft DAB	270462	1
Gehe	Zinc paste, soft DAB	1B110112	1
GEHE	Hydrogen peroxide sol...	H216091/1	1
GEHE	Hydrogen peroxide sol...		15
Gehe 06.08.2013	Hydrogen peroxide sol...	6312065	1
Gehe 06.09.13	Rosemary leaf tinctur...	11323611	1
gehe 11.06.2013	Zinc paste DAB	276141	2
gehe 11.06.2013	Zinc paste, soft DAB	12D10-T03	1
Gehe 12.04.2013	Lime water	12177901	1
Gehe 16.04.2013	Coconut oil, hardened	13010803	1
Gehe 16.05.2013	Comfrey root tincture...	12054016	2
Gehe 26.09.13	Wormwood tincture (70...	12237913	2
Gehe/Fischar/11.12.2012	Hydrogen peroxide sol...	6312105	8
GEHE/Fischar/12.06.2012	Hydrogen peroxide sol...	6312012	3
Gehe/Fischar11.07.2012	Hydrogen peroxide sol...	6312051	9
Gehe/Maros	Opium tincture, stand...	11006112	17
Gehrlicher	Pansy mother tincture	RU1775	1
Gehrlicher/Gehrlicher	Burnet root mother ti...	U3600	1
Großhandel	Arnica tincture (70%)...	1204046-01	4
Großhandel	Arnica tincture (70%)...	1211003-01	5
Großhandel	Coal tar solution	111031-01	1
Großhandel	Sage tincture (70%) 1...	11220912	1
Großhandel	Sage tincture (70%) 1...	12196612	5
Großhandel	Vitamin A-palmitate (...)	13057007	4
Großhandel	Zinc oxide lotion, aq...	13087501	13
Großhandel	Zinc paste DAB	12147903	2
Großhandel	Zinc paste DAB	13043502	4
H.Lamotte/Noweda	Castor oil, hydrated	123384	1
Hanosan	Strophanthus seed mot...	316899	1
Hanosan	Strophanthus seed mot...	348059	1
Hedinger	Formic acid 85%	383338	1
Hedinger	Formic acid 85%	170912	1
Hedinger	Hydrogen peroxide sol...	301208	1
Hedinger	Hydrogen peroxide sol...	301205	1
Hedinger	Hydrogen peroxide sol...	301212	1
Hedinger	Hydrogen peroxide sol...	301207	1
Hedinger	Hydrogen peroxide sol...	301211	1
Hedinger	Hydrogen peroxide sol...	301210	1
Hedinger	Hydrogen peroxide sol...	61011b	2
Hedinger	Hydrogen peroxide sol...	61008b	1
Hedinger	Hydrogen peroxide sol...	610108	1
Hedinger	Hydrogen peroxide sol...	61010b	2
Hedinger	Hydrogen peroxide sol...	8111302	1
Hedinger	Hydrogen peroxide sol...	210913H202	1
Hedinger	Hydrogen peroxide sol...	61006B	2
Hedinger	Hydrogen peroxide sol...	61003B	2
Hedinger	Hydrogen peroxide sol...	16092013A	1
Hedinger	Hydrogen peroxide sol...	61005B	3
Hedinger	Hydrogen peroxide sol...	610128	1
Hedinger	Witch hazel bark water	20153B	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
HEDINGER	Hydrogen peroxide sol...	601008B	1
Hedinger / Gehe	Hydrogen peroxide sol...	64044B	1
Hedinger / Gehe	Hydrogen peroxide sol...	(061b)061b	2
Hedinger / Sanacorp	Sodium hypochlorite s...	46386	1
Hedinger Jenne	Hydrogen peroxide sol...	61011B	28
Hedinger Jenne	Hydrogen peroxide sol...	61009B	24
Hedinger/ Gehe	Sodium hypochlorite s...	46375	1
Hedinger/ Phönix	Sodium hypochlorite s...	46392	1
Hedinger/AHD	Hydrogen peroxide sol...	61010B	1
Hedinger/AHD	Hydrogen peroxide sol...	61012B	8
Hedinger/AHD	Hydrogen peroxide sol...	61010B	26
Hedinger/Anzag	Hydrogen peroxide sol...	4413A-01720	1
Hedinger/Anzag 21.07...	Hydrogen peroxide sol...	61008B	21
Henry Lamotte	Zinc oxide lotion, aq...	60480	1
Henry Lamotte GmbH/Ho...	Zinc oxide lotion, aq...	5289300	3
henry lamotte oils gmbh	Zinc oxide lotion, aq...	6000083002	1
henry lamotte oils gmbh	Zinc oxide lotion, aq...	90616-1	1
Henry Lamotte Oils GmbH	Zinc oxide lotion, aq...	600320602	1
Henry Lamotte Oils GmbH	Zinc paste, soft DAB	5316600	8
Henry Lamotte Oils Gm...	Zinc paste, soft DAB	5375700	4
Henry Lamotte Oils, A...	Zinc oxide lotion, aq...	5375600	12
Henry Lamotte/ 24.04...	Zinc oxide lotion, aq...	1221401	2
HenryLamotte/Sanacorp	St. John's wort oil	5308402-3457312120	1
Herbeta Arzneimittel	Hydrogen peroxide sol...	15061202	1
Herstellung	Hydrogen peroxide sol...	6313072	1
Hetterich/Teofarma	Myrrh tincture (90%) ...	139101	1
Holdermann/Lamotte	Castor oil, hydrated	5280301	2
Hubertus Apo	Hydrogen peroxide sol...	240513	1
Hubertus Apo	Hydrogen peroxide sol...	100913	1
Hubertus Apotheke	Hydrogen peroxide sol...	5082014	1
Hubertus Apotheke	Hydrogen peroxide sol...	130114	1
Hustadt Apotheke	Hydroxyethylcellulose...	12128101	1
Hustadt Apotheke	Ratanhia tincture (70...	12255401	7
ichthyol	Clobetasol 0.5% Corde...	5112M-03390	1
ichthyol	Clobetasol 0.5% Corde...	1713M-03390	1
Ichthyol	Betamethason-V 1.22% ...	11D01	1
Ichthyol	Betamethason-V 1.22% ...	11/2015	1
Ichthyol	Betamethason-V 1.22% ...	110001	1
Ichthyol	Clobetasol 0.5% Corde...	11D003	15
Ichthyol	Clobetasol 0.5% Corde...	12d012	1
Ichthyol	Clobetasol 0.5% Corde...	10D024	1
Ichthyol	Clobetasol 0.5% Corde...	9D035	1
Ichthyol	Methoxsalen 0.006% Co...	130004	2
Ichthyol	Methoxsalen 0.006% Co...	14D001	1
Ichthyol- Gesellschaft	Methoxsalen 0.006% Co...	13D013	7
Ichthyol Gesellschaft...	Clobetasol 0.5% Corde...	11D034	25
ICHTHYOL/ahp	Methoxsalen 0.006% Co...	13d004	1
Ichthyol/Gehe	Methoxsalen 0.006% Co...	13D004	32
Ichthyol/Jenne	Betamethason-V 1.22% ...	12D015	26
Ichthyol/Jenne	Betamethason-V 1.22% ...	12D015	2
Ichthyol/Noweda	Betamethason-V 1.22% ...	120015	1
Ichthyol-Gesellschaft...	Betamethason-V 1.22% ...	11D001	58
Ichthyol-Gesellschaft...	Gentamicin 1% Cordes®...	9D002	1
Ichtyol	Clobetasol 0.5% Corde...	12D012	13
Infectopham/Gehe	Hydrophobic triclosan...	S021201.	1
infectopharm	Hydrophobic triclosan...	S041302.1	2
infectopharm	Hydrophobic triclosan...	s101209.1	1
Infectopharm	Hydrophobic triclosan...	S031101.1	1
Infectopharm	Hydrophobic triclosan...	S091203.1	8
Infectopharm	Hydrophobic triclosan...	20160712-1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Infectopharm	Hydrophobic triclosan...	S101208.1	14
Infectopharm	Hydrophobic triclosan...	S031502.1	1
Infectopharm	Hydrophobic triclosan...	S0212041	1
Infectopharm	Hydrophobic triclosan...	S041301.1	7
Infectopharm	Hydrophobic triclosan...	S091109.1	2
Infectopharm	Hydrophobic triclosan...	S101207.1	10
Infectopharm	Hydrophobic triclosan...	S04131011	1
Infectopharm	Hydrophobic triclosan...	s1012091	1
InfectoPharm	Hydrophobic triclosan...	S091103.1	3
InfectoPharm / 04.09...	Hydrophobic triclosan...	S091104.1	1
InfectoPharm / Gehe	Hydrophobic triclosan...	S0413041	1
Infectopharm, Noweda	Hydrophobic triclosan...	S021203.1	4
Infectopharm, Noweda	Hydrophobic triclosan...	S0231204.1	1
Infectopharm,AHC	Hydrophobic triclosan...	S101304.1	1
Infectopharm/AHD	Hydrophobic triclosan...	S1012209.1	1
InfectoPharm/GEHE	Hydrophobic triclosan...	S.091203.1	1
infectopharm/Noweda	Hydrophobic triclosan...	S041301.1	5
Infectopharm/Noweda	Hydrophobic triclosan...	S0413011	1
InfectoPharm/Phoenix	Hydrophobic triclosan...	S021201.1	4
InfectoPharm/Phönix	Hydrophobic triclosan...	1101410005	1
Jenne	Myrrh tincture (90%) ...	12197205	5
Kehr	Zinc ointment DAB	272144	1
Kehr	Zinc ointment DAB	274127	2
Kehr	Zinc oxide lotion, aq...	12065202	2
Klenk	Calendula oil, fat	12308503	2
Krieger	Coal tar solution	11319212	3
labor	Hydrogen peroxide sol...	7571	1
Labor	Hydrogen peroxide sol...	7408	1
Labor	Hydrogen peroxide sol...	9020	1
Labor	Hydrogen peroxide sol...	9266	1
Labor	Hydrogen peroxide sol...	9758	1
Labor	Hydrogen peroxide sol...	8762	1
Labor	Hydrogen peroxide sol...	9848	1
Labor	Hydrogen peroxide sol...	8698	1
Labor	Hydrogen peroxide sol...	9634	1
Labor	Hydrogen peroxide sol...	7991	1
Labor	Hydrogen peroxide sol...	7926	1
Labor	Hydrogen peroxide sol...	7143	1
Labor	Hydrogen peroxide sol...	7799	1
Labor	Hydrogen peroxide sol...	7798	1
Labor	Hydrogen peroxide sol...	9957	1
Labor	Hydrogen peroxide sol...	9757	1
lamotte	Castor oil, hydrated	5188401	2
lamotte	Zinc oxide lotion, aq...	6003206002	2
Lamotte	Arnica tincture (70%)...	5139400	1
Lamotte	St. John's wort oil	5346702	14
Lamotte	Zinc oxide lotion, aq...	5651500	1
Lamotte	Zinc oxide lotion, aq...	5099300	1
Lamotte	Zinc oxide lotion, aq...	5735500	3
Lamotte	Zinc oxide lotion, aq...	5163600	2
Lamotte	Zinc oxide lotion, aq...	60460	1
Lamotte	Zinc oxide lotion, aq...	5216200	1
Lamotte	Zinc paste, soft DAB	5216300	4
Lamotte / Gehe	Castor oil, hydrated	5138701	3
Lamotte / Gehe	Zinc oxide lotion, aq...	5517400	12
Lamotte / Gehe	Zinc oxide lotion, aq...	5321300	10
Lamotte / Sanacorp	St. John's wort oil	517601	1
Lamotte Sanacorp	St. John's wort oil	5308402	4
Lamotte/ Finze	Zinc oxide lotion, aq...	5563300	2
lamotte/Noweda	Zinc paste, soft DAB	5447900	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Leyhs/Noweda	Methylrosanilinium ch...	124141	1
Maros	Angelica root tinctur...	1028713	2
Maros	Angelica root tinctur...	8030113	1
Maros	Bitter candytuft fres...	5408	1
Maros	Bitter candytuft fres...	5433	1
Maros	Marian thistle fruit ...	U3762	2
Maros	Matricaria liquid ext...	6990	3
Maros	Melissa leaf liquid e...	6399	3
Maros	Opium tincture, stand...	6006713	32
Maros	Opium tincture, stand...	4005612	4
Maros	Opium tincture, stand...	3009616	1
Maros	Opium tincture, stand...	2009416	2
Maros	Opium tincture, stand...	9006012	16
Maros	Opium tincture, stand...	7008915	1
Maros	Opium tincture, stand...	8005812	7
Maros	Opium tincture, stand...	13071303	1
Maros	Opium tincture, stand...	10009115	2
Maros	Opium tincture, stand...	2004511	1
Maros	Opium tincture, stand...	1006313	14
Maros	Opium tincture, stand...	3007514	2
Maros	Opium tincture, stand...	4009716	1
Maros	Opium tincture, stand...	8009015	2
Maros	Opium tincture, stand...	15031302	1
Maros	Peppermint leaf mothe...	RU2432	2
Maros	Peppermint tincture (...)	RU2432	1
Maros	St. John's wort oil	9758	1
Maros	St. John's wort oil	97583x	1
Maros	Thyme liquid extract	1156	1
MAROS	Opium tincture, stand...	3005512	1
Maros Arznei GmbH	Opium tincture, stand...	4006513	34
Maros Arznei GmbH	Opium tincture, stand...	8006913	27
Maros Arznei GmbH	Opium tincture, stand...	11007113	8
MAROS Arznei GmbH	Opium tincture, stand...	7006813	17
Maros Arznei GmbH/Anzag	Opium tincture, stand...	2005412	1
Maros/direkt/07.05.13	Strophanthus seed mot...	4029013	1
Maros/direkt/15.01.14	Strophanthus seed mot...	1031014	3
Maros/direkt/16.08.13	Strophanthus seed mot...	8029813	1
Maros/direkt/21.10.13	Strophanthus seed mot...	10030413	2
Maros/direkt/28.02.14	Strophanthus seed mot...	2031214	5
Maros/Maros	Cistus tincture 1:5	11027812	1
Maros/Maros	Strophanthus seed mot...	6029813	2
Maros/Phoenix	Opium tincture, stand...	60066713	1
müller	Propolis tincture (96%)	131206	1
Nachprüfung	Chinese cinnamon oil		1
Nachprüfung	Gentian tincture (70%...	411E02442	1
Nachprüfung	Peru balsam, original	12020602	1
Nachprüfung	Standardised cayenne ...	2413M-03232	1
noweda	Witch hazel bark water	12100302	1
Noweda	Arnica tincture (70%)...	1106037-01	1
Noweda	Betamethason-V 1.22% ...	11D035	104
Noweda	Castor oil, hydrated	12114603	4
Noweda	Castor oil, hydrated	268877	2
Noweda	Coal tar solution	11319204	10
Noweda	Gentamicin 1% Cordes® ...	12D010	28
Noweda	Gentamicin 1% Cordes® ...	10D026	35
Noweda	Methoxsalen 0.006% Co...	12D011	9
Noweda	Methoxsalen 0.006% Co...	11D006	1
Noweda	Myrrh tincture (90%) ...	1219104	1
Noweda	Opium tincture, stand...	1106112	2
Noweda	Orange flower oil, ar...	10267805	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Ratanhia tincture (70...	12255406	1
Noweda	Thyme liquid extract	10198605	1
Noweda	Triamcinolone acetoni...	12080701	1
Noweda	Triamcinolone acetoni...	11359001	1
Noweda	Valerian tincture (70...	F1300052	1
Noweda	Valerian tincture (70...	F1300146	1
Noweda	Vitamin A-palmitate (...)	13F26-N08	5
Noweda	Vitamin A-palmitate (...)	13C19-N12	1
Noweda	Witch hazel bark water	11062805	1
Noweda	Witch hazel leaf liqu...	12095405	1
Noweda	Zinc ointment DAB	11271001	3
Noweda	Zinc oxide lotion, aq...	11329903	3
Noweda	Zinc oxide lotion, aq...	11281401	2
Noweda	Zinc oxide lotion, aq...	1211701	1
Noweda	Zinc paste DAB	276606	1
Noweda	Zinc paste, soft DAB	11338501	3
Noweda	Zinc paste, soft DAB	1138601	1
Noweda/apomix	Clobatesol propionate...	13A131113	1
Noweda/Fischar	Hydrogen peroxide sol...	6313024	1
Noweda/Fischar	Hydrogen peroxide sol...	6313052	1
Noweda/Fischar	Hydrogen peroxide sol...	6313075	2
Noweda/Lamotte	Zinc oxide lotion, aq...	5447800	9
Otto Fischar	Hydrogen peroxide sol...	6313101	1
Otto Fischar / Noweda	Hydrogen peroxide sol...	6314011	1
Otto Fischar gmbh	Hydrogen peroxide sol...	6315031	1
Otto Fischar gmbh	Hydrogen peroxide sol...	6314113	1
Otto Fischar GmbH	Hydrogen peroxide sol...	6312021	3
Otto Fischar GmbH	Rose water	5354081	1
Otto Fischar GmbH & C...	Lactic acid 80% L(+)	20020271	1
Otto Fischar GmbH& Co...	Hydrogen peroxide sol...	6313055	5
Park Apotheke	Matricaria flower tin...	nachprüfung	1
pharmachem	Zinc oxide oil DAC	90512	1
pharmachem	Zinc oxide paste, sof...	211112	1
Pharmachem	Oxytetracycline ointm...	230113	1
Pharmachem	Zinc oxide cream 25% ...	110112	1
Pharmachem	Zinc oxide oil SR 90	150413	1
Pharmachem	Zinc oxide oil SR 90	160413	1
Pharmachem	Zinc oxide oil SR 90	30712	1
Pharmachem	Zinc oxide oil SR 90	20712	1
Pharmachem	Zinc oxide oil SR 90	260412	1
Pharmachem	Zinc oxide paste, sof...	70512	1
PharmaChem	Zinc oxide cream 25% ...	150513	1
Pharmachem	Oxytetracycline ointm...	140212	2
Pharmachem Anzag	Zinc oxide oil SR 90	140513	1
Pharmachem Anzag	Zinc paste, soft DAB	191112	1
Pharmachem GmbH	Zinc oxide oil SR 90	61112	1
PharmaChem/Noweda	Zinc oxide oil DAC	80512	1
Pharmachem/Phönix	Zinc oxide oil SR 90	160810	1
pharmachen	Zinc oxide oil SR 90	30045100	1
Pharmachern/AHD	Zinc oxide oil SR 90	160122	1
Phoenix	Castor oil, hydrated	12071107	1
Phoenix	Castor oil, hydrated	11325205	2
Phoenix	Zinc oxide lotion, aq...	4010501	1
Phoenix	Zinc paste, soft DAB	281120	1
Phoenix (Pharmachem P...	Zinc oxide paste, sof...	191112	1
Phoenix (PKH Halle)	Methylrosanilinium ch...	4A200712	1
Phoenix (PKH Halle)	Metronidazole 10% (no...	11B251012	1
Phoenix 19.12.12	Vitamin A-palmitate (...)	12J16-N03	2
Phoenix am 19.09.2012	Coconut oil, hardened	11379505	1
Phoenix, 17.09.2012	St. John's wort oil	12076905	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix/Apomix	Clobatesol propionate...	6A070513	1
phönix	Camphor oil 10%	7/AKD57	2
Phönix	Calendula tincture (7...	13142803	10
Phönix	Calendula tincture (7...	13142811	9
Phönix	Cayenne pepper tinctu...		1
Phönix	Coal tar solution	12105-n02	1
Phönix	Coal tar solution	11B14-N04	1
Phönix	Hydrogen peroxide sol...	6313073	2
Phönix	Hydroxyethylcellulose...	13286702	1
Phönix	Matricaria liquid ext...	12360301	1
Phönix	Opium tincture, stand...	2006413	28
Phönix	Sage leaf liquid extr...	11372206	4
Phönix	Siam benzoin tincture...	12228107	1
Phönix	St. John's wort oil	5803902	2
Phönix	St. John's wort oil	12293904	4
Phönix	Theriac	12067210	1
Phönix	Thyme liquid extract	12197005	3
Phönix	Thyme tincture (70%) ...	10379114	1
Phönix	Triamcinolone acetoni...	13029801	2
Phönix	Zinc oxide lotion, aq...	11329801	3
Phönix	Zinc oxide oil DAC	286040	1
Phönix	Zinc oxide oil DAC	10010212	1
Phönix	Zinc paste, soft DAB	15005101	1
Phönix	Castor oil, hydrated	288639	1
Phönix 11.11.2014	Food colorant green ...	13020106	1
Phönix Gotha	Castor oil, hydrated	12H29-N04	2
Phönix Gotha	Castor oil, hydrated	12H29-N04	1
Phönix Gotha	Gentian tincture (70%...	10334601	1
Phönix Gotha	Hydrophilic metronida...	4B190413	6
Phönix Gotha	Hydrophobic triclosan...	S091201.1	1
Phönix Gotha	Hydrophobic triclosan...	561069	1
Phönix Gotha	Hydroxyethylcellulose...	11341802	1
Phönix Gotha	Hydroxyethylcellulose...	12320501	8
Phönix Gotha	Hydroxyethylcellulose...	12244003	5
Phönix Gotha	Lactic acid 80% L(+)	2002071	1
Phönix Gotha	Peanut oil, hydrated	12230701	2
Phönix Gotha	Peru balsam, original	10C30-N07	1
Phönix Gotha	Valerian tincture (70...	11254302	1
Phönix Gotha	Vitamin A-palmitate (...)	8900521	1
Phönix Gotha	Wormwood tincture (70...	10010909	1
Phönix Gotha	Zinc oxide cream 25% ...	276963	1
Phönix Gotha	Zinc oxide lotion, aq...	259142	1
Phönix Gotha	Zinc paste, soft DAB	113381	1
Phönix Maros	Castor oil, hydrated	285801	1
Phönix Maros	Opium tincture, stand...	11008214	2
Phönix Maros	Opium tincture, stand...	7007914	2
Phönix, 16.10.2012	Coal tar solution	12E08-N01	5
phönix, 21.11.2012	Rose water	12C27-N05	2
phönix, 23.11.2012	Hydroxyethylcellulose...	24B070912	2
Phönix, 30.10.2012	Valerian tincture (70...	12113005	2
Phönix/ henry lamotte	St. John's wort oil	5442801	1
Phönix/ henry lamotte	St. John's wort oil	5171601	7
Phönix/ Maros	Opium tincture, stand...	5006613	30
Phönix/Eurho	Castor oil, hydrated	1212008-02	1
Phönix/euro	Castor oil, hydrated	1504006-02	1
Phönix/Henry Lamotte	St. John's wort oil	5540603	2
Phönix/HenryLaMotte	Zinc paste, soft DAB	56119000	1
Phönix/InfectoPharm	Hydrophobic triclosan...	S021204.1	7
Phönix/InfectoPharm	Hydrophobic triclosan...	S021205.1	5
PKH / Gehe	Hydrogen peroxide sol...	13A02-M09	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PKH / Gehe	Hydrogen peroxide sol. ...	12I06-M07	2
PKH / Gehe	Hydrogen peroxide sol. ...	12I04-M09	1
PKH / Gehe	Hydrogen peroxide sol. ...	12D27-M09	1
PKH / Gehe	Hydrogen peroxide sol. ...	12D26-M05	1
PKH GbmH	Betamethasone valerat. ...	3B280213	2
PKH GbmH	Clobatesol propionate. ...	2A280213	1
PKH GmbH	Hydrogen peroxide sol. ...	487A161012	1
PKH GmbH	Zinc oxide cream 25% ...	2B080113	1
PKH GmbH/Anzag	Hydrogen peroxide sol. ...	485A121012	1
PKH Halle	Searing paste SR	17B211112	1
PKH Halle / Phoenix	Metronidazole 10% (no. ...)	8B211013	1
PKH Halle / Phoenix	Searing paste SR	14B240913	4
PKH Halle / Phoenix	Searing paste SR	5B110413	2
PKH Halle/ Alliance H.	Hydroxyethylcellulose. ...	18B090813	4
Post-Apotheke	Food colorant green ...	13020101	1
primavera	St. John's wort oil	12106131	1
Primavera	Aloe vera oil	8601131	1
Primavera	St. John's wort oil	25508131	1
Purren Apotheke	St. John's wort oil	11205311	1
Resena/Anzag	Carotene (ca. 4m IU)	1204-3534	1
Retest	Coal tar solution	11D01-N04	1
s.Prüfprotokoll	Peru balsam, original	s.Prüfprotokoll	1
sanacorp	Opium tincture, stand. ...	9005912	5
Sanacorp	Benzalkonium chloride	10190002	1
Sanacorp	Camphor oil 10%	7/AKD4	1
Sanacorp	Hydrophobic triclosan. ...	S031105.1	1
Sanacorp	Hydroxyethylcellulose. ...	12127901	2
Sanacorp	Hydroxyethylcellulose. ...	12309303	9
Sanacorp	Hydroxyethylcellulose. ...	13190501	8
Sanacorp	Hydroxyethylcellulose. ...	12219401	6
Sanacorp	Opium tincture, stand. ...	60006713	1
Sanacorp	Orange flower water	AKD34	1
Sanacorp	St. John's wort oil	13252505	4
Sanacorp	Witch hazel leaf liqu. ...	13112003	1
Sanacorp	Zinc oxide lotion, aq. ...	AKD33	1
Sanacorp	Zinc paste DAB	13043401	21
Sanacorp	Zinc paste DAB	13201903	1
Sanacorp	Zinc paste DAB	13345101	9
Sanacorp	Zinc paste, soft DAB	274978	3
Sanacorp	Zinc paste, soft DAB	274979	1
Sanacorp, 1303.12 EK: ...	Hydrophobic triclosan. ...	S091102.1	1
Sanacorp/ apomix	Zinc oxide oil DAC	26B190911	1
Sanacorp/ PKH Halle	Hydroxyethylcellulose. ...	10B200412	1
Sanacorp/ PKH Halle	Hydroxyethylcellulose. ...	8B050412	1
Sanacorp/Ichthyol	Methoxsalen 0.006% Co. ...	13D013	1
Sigmar Aldrich	Dichloroacetic acid 10%	STBC3539V	1
Sino Phyto	Hydrogen peroxide sol. ...	8517	1
Spinnrad/Phoenix	Witch hazel bark water	413310062	1
stumpf	Hydrogen peroxide sol. ...	280A060712	1
Taoasis	Aloe vera oil	130501-112371	1
Taoasis	Aloe vera oil	130237-112020BAG94	1
Taoasis	Aloe vera oil	130237-111659	1
Taoasis	Aloe vera oil	111208-105327	1
Taoasis	Aloe vera oil	120247-107323	1
Taoasis	Aloe vera oil	120247-106498	1
Taoasis	Aloe vera oil	110137-96542	1
Taoasis	Arnica oil, fat	7360-112511	1
Taoasis	Arnica oil, fat	7360-111459	1
Taoasis	Arnica oil, fat	7360-111366	1
Taoasis	Arnica oil, fat	7525-113064	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Taoasis	Calendula oil, fat	119943-110316	1
Taoasis	Calendula oil, fat	119426-108734	1
Taoasis	Siam benzoin	541389-111081	1
Taoasis	St. John's wort oil	145646-103095	1
Taoasis	St. John's wort oil	123212-111876	1
Taoasis	St. John's wort oil	123212-110835	1
Taoasis	St. John's wort oil	12312-109683	1
Taoasis GmbH	Aloe vera oil	130237	1
Taoasis GmbH	St. John's wort oil	150822-107470	1
Taoasis	Macadamia oil, refined	4976601-92844	1
Teofarma	Arnica tincture (70%)...	144301	2
Wasserfuhr	Zinc paste DAB	12169401	1

Appendix D: Requirements of validation

In order to ensure adherence to the safe scientific status, the individual methods for manufacturing and testing must be validated under certain circumstances (compare § 34 para. 1 no. 3, § 35 para. 1 no. 4 and para. 4 sentence 1 no. 2 b, para. 6 sentence 3 *ApoBetrO* [Pharmacies Rules and Regulations]). The *ApoBetrO* [Pharmacies Rules and Regulations] incorporates a legal definition in § 1 a para. 16 (quotation translated):

“Validation is the provision of documented proof which with a high degree of safety documents that, via a specific process or standard work process, a medicinal product is manufactured and tested, which is in accordance with previously determined quality features.”

Validation documentation can be used to prove that methods or devices which are not described in the Pharmacopoeia within the meaning of § 6 para. 1 sentence 3 *ApBetrO* [Pharmacies Rules and Regulations] achieve the same results as those in the Pharmacopoeia. On the other hand, with the requirements of the demanded validation it must be observed whether the respective testing method is already incorporated in the Pharmacopoeia.

NIR spectroscopy as a general testing method need not be validated in accordance with the express ruling in the *Ph. Eur. Section 1.1* [3], as it is already described in *Section 2.2.40* of the *Ph. Eur.* as an area of application for the identification of raw materials.

However, a special validation requirement exists for the reference database. This requirement is met with the existing document. Further requirements or rules as to how this proof must be furnished do not exist. It is required that the processes guarantee the same results as the methods and devices in the Pharmacopoeia [17].

Carrying out identity tests with *Apo-Ident* is therefore also possible if the NIR spectroscopy process is not required in the Pharmacopoeia monograph of the substance for identity testing. All NIR analyses with *Apo-Ident* prove several, often all molecule groups and are therefore comparable with a series of individual, targeted chemical proofs [4]. Therefore, the identity proof with *Apo-Ident* replaces the monograph test series (with two or more test combinations).

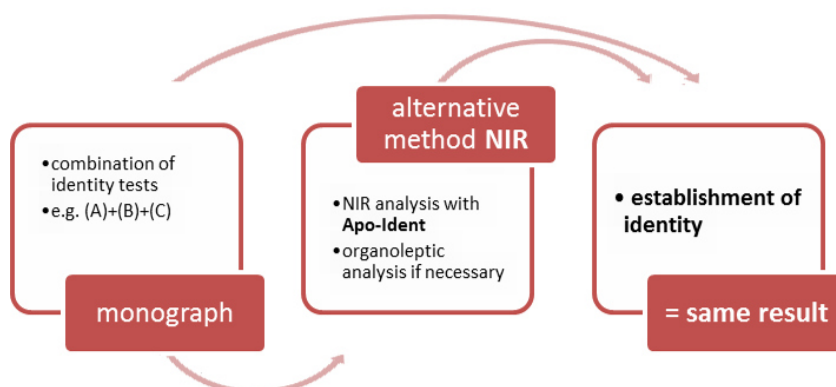


Figure 2: The combination of tests of the monograph is replaced by the alternative method NIR spectroscopy using *Apo-Ident*. This is permissible because both test procedures result in the establishment of the identity of the raw material.

With this validation documentation, proof is furnished that identical results are achieved with *Apo-Ident* and Pharmacopoeia methods, i.e. confirmation of the identity of the raw material [2].

Appendix E: Conformity of Apo-Ident with the European pharmacopeia

According to *Ph. Eur. Section 2.2.40*, NIR spectroscopy is basically suitable for: “Identification of agents, excipients, dosage forms, intermediate manufacturing products, chemical raw materials and packaging materials” ([3], quotation translated).

The fact that *Apo-Ident* meets the further criteria of the European Pharmacopoeia under the headings in *Section 2.2.40*

- Apparatus
- Measurement methods
- Sample preparation and presentation
- Testing the functionality of the instrument
- Identification and characterisation (qualitative analysis)
- Quantitative analysis
- Ongoing model evaluation
- Transfer of databases
- Data storage

can be proven based on the *HiperScan GmbH* documentation of “Meeting 2.2.40 *Ph. Eur.* by *Apo-Ident*” [4].

Appendix F: Addition to the second-stage models

As described in the general part of the *Apo-Ident* validation documentation, in the section [Model creation procedure and validation runs](#), under 2. *Generating the chemometric models (calibration) in f)*, it is possible to create a further chemometric model for certain substances (second-stage model) and to undertake the assessment in several steps.

Substances in second-stage models

In the *Apo-Ident Update 2019-01*, a distinction is made between the following substances in the listed second-stage models:

Submodel 1

Bacitracin
Betamethasone, micronized
Betamethasone valerate
Budesonide, micronized
Capsaicin
Dexamethasone
Erythromycin
Heparin sodium
Hydrocortisone butyrate
Norethisterone acetate
Prednicarbate, micronized
Prednisolone, micronized
Prednisone
Triamcinolone, micronized
Triamcinolone acetonide

Submodel 2

Beclometasone dipropionate, anhydrous
Betamethasone dipropionate, micronized
Clobetasol propionate
Diphenylcyclopropenone
Gentamicin sulfate
Sodium benzoate
Sodium citrate
Oxybutynin hydrochloride

Submodel 3

Alfatradiol
Estriol
Pregnenolone
Progesterone, micronized
Quinine hydrochloride
Quinine sulfate dihydrate
Spironolactone
Testosterone propionate

Liquid submodel

2-ethylhexyl laurate
Base cream DAC

Lanette® ointment (conserved)
Non-ionic hydrophil ointment
Squalane
Wool alcohols cream DAB/SR

References

- [1] ABDA – BUNDESVEREINIGUNG DEUTSCHER APOTHEKERVERBÄNDE: Verordnung über den Betrieb von Apotheken (Apothekenbetriebsordnung – ApBetrO), 2012
- [2] REIMANN, B. ; REGIERUNGSPRÄSIDIUM DARMSTADT: Hinweise zur ordnungsgemäßen Prüfung von Arzneimitteln und Ausgangsstoffen (§§ 6 und 11 *ApBetrO*), 2007
- [3] *Europäisches Arzneibuch, Grundwerk 2014 einschließlich 1. bis 4. Nachtrag.* 8. Ausgabe. Deutscher Apotheker Verlag (978-3-7692-6508-8)
- [4] HIPERSCAN GMBH: Erfüllung von 2.2.40 Ph. Eur. durch Apo-Ident, 2013
- [5] ARBEITSGEMEINSCHAFT DER PHARMAZIERÄTE DEUTSCHLANDS (APD): Resolution in Verbindung mit §§ 6 und 11 *ApBetrO* Verwendung eines Nah-Infrarot-Spektrometers (NIR) zur Identitätsprüfung, 16. 10. 2013, DAZ 48, November 2013
- [6] ARBEITSGEMEINSCHAFT DER PHARMAZIERÄTE DEUTSCHLANDS (APD): Resolution 2014, Arbeitsgemeinschaft der Pharmazieräte Deutschlands (APD), Oktober 2014
- [7] *DAC/NRF.* Govi-Verlag (978-3-7741-0044-2)
- [8] KESSLER, W.: *Multivariate Datenanalyse.* WILEY-VCH Verlag, 2007 (978-3-527-31262-7)
- [9] NÆS, T. ; ISAKSSON, T. ; FEARN, T. ; DAVIES, T.: *Multivariate Calibration and Classification.* NIR Publications, 2002 (978 0 9528666 2 6)
- [10] HANLEY, J. A. ; LIPPMAN-HAND, A.: If nothing goes wrong, is everything all right? In: *Journal of the American Medical Association* 249 (1983), S. 1743–1745
- [11] JOVANOVIĆ, B. D. ; LEVY, P. S.: A Look at the Rule of Three. In: *The American Statistician* 51 (1997), S. 137–139
- [12] BRONSTEIN, I. N. ; SEMENDJAJEW, K. A. ; MUSIOL, G. ; MÜHLIG, H.: *Taschenbuch der Mathematik.* 5. überarbeitete und erweiterte Auflage. Verlag Harri Deutsch, 2000 (3-8171-2015-2)
- [13] MAHALANOBIS, P.: On the generalized distance in statistics. In: *Proc. Nat. Inst. Sci. India (Calcutta)* 2 (1936), S. 49–55
- [14] YAMBOR, B. ; DRAPER, W. ; BEVERIDGE, R.: Analyzing PCA-based face recognition algorithms: Eigenvector selection and distance measures. In: *Second Workshop Empirical Evaluation in Computer Vision* (2000)
- [15] HIPERSCAN GMBH: Identifikationsmethodik Apo-Ident, 2012
- [16] HIPERSCAN GMBH: Datenvorbehandlung des Identifikationssystems Apo-Ident, 2012
- [17] CYRAN ; RÖTTA: Apothekenbetriebsordnung, Kommentar § 6 Rn. 10, 2010

Index

- 1,3-Butandiolum, *see* [1,3-butanediol](#)
 1,3-butanediol, [1223](#)
 2-Ethylhexylislauromyristas, *see* [2-ethylhexyl laurate](#)
 2-ethylhexyl laurate, [1227](#)
- Acetic acid 30%, [21](#)
 Acetic acid 99%, [25](#)
 Acetone, [31](#)
 Acetonum, *see* [Acetone](#)
 Acidum aceticum 30 %, *see* [Acetic acid 30%](#)
 Acidum aceticum dilut. 30%, *see* [Acetic acid 30%](#)
 Acidum aceticum glaciale, *see* [Acetic acid 99%](#)
 Acidum glycolicum 70%, *see* [Glycolic acid 70%](#)
 Acidum L (+)-lacticum 20%, *see* [Dextrorotary lactic acid 20%](#)
 Acidum lacticum, *see* [Lactic acid](#)
 Acidum oleicum, *see* [Oleic acid](#)
 Acidum oleinicum, *see* [Oleic acid](#)
 Acidum salicylicum cum vaselino alba 50%, *see* [Salicylic acid 50% \(paraffin\)](#)
 Acorus oil, [41](#)
 Adeps lanae, *see* [Wool fat](#)
 Adeps lanae anhydricus, *see* [Wool fat](#)
 Adeps solidus 378, *see* [Softisan® 378](#)
 Adeps solidus 649, *see* [Hard fat 649](#)
 Aether aceticus, *see* [Ethyl acetate](#)
 Aether petrolei, *see* [Petroleum ether](#)
 Alcohol isopropylicus 100%, *see* [Isopropyl alcohol 100%](#)
 Alcohol isopropylicus 70%, *see* [Isopropyl alcohol 70%](#)
 Alkohol benzylicus, *see* [Benzyl alcohol](#)
 All-rac-alpha-tocopheryl acetate, [45](#)
 Almond oil, refined, *see* [Triglycerides](#)
 alpha-Tocopherolum, *see* [D-alpha tocopherol](#)
 Althae sirupus (konserviert), *see* [Syrups](#)
 Aluminium acetotartrate solution, *see* [Water and aqueous solutions](#)
 Amygdalae oleum raffinatum, *see* [Triglycerides](#)
 Angelica archangelica oleum, *see* [Angelica root oil](#)
 Angelica root oil, [51](#)
 Anionic hydrophilic cream DAB, *see* [Hydrophilic ointments](#)
 Anionic hydrophilic cream SR (NRF p.27), *see* [Hydrophilic ointments](#)
 Anionic hydrophilic cream with propylene glycol, *see* [Hydrophilic ointments](#)
 Anionic hydrophilic cream without preservatives, *see* [Hydrophilic ointments](#)
 Anise oil, *see* [Anise oil / Star aniseed oil](#)
 Anise oil / Star aniseed oil, [55](#)
 Apricot kernel oil, refined, *see* [Triglycerides](#)
 Aqua, *see* [Water and aqueous solutions](#)
 Aqueous salicylic acid 50% (non-ionic, hydrophilic cream SR), [61](#)
 Arachidis oleum raffinatum, *see* [Triglycerides](#)
 Avocado oil, refined, *see* [Triglycerides](#)
 Avocado oleum raffinatum, *see* [Triglycerides](#)
- Base cream DAC, [67](#)
 Base cream, hydrophobic DAC, *see* [Hydrophilic ointments](#)
 Basis Cordes® RK, [121](#)
 Benzalkonii chloridi solutio 50%, *see* [Benzalkonium chloride solution 50%](#)
 Benzalkonium chloride solution 50%, [131](#)
 Benzinum (petrolei), *see* [Petroleum ether](#)
 Benzyl alcohol, [137](#)
 Benzyl benzoate, [141](#)
 Benzyl nicotinate, [145](#)
 Benzylis nicotinas, *see* [Benzyl nicotinate](#)
 Benzylum benzoicum, *see* [Benzyl benzoate](#)
 Benzylum nicotinicum, *see* [Benzyl nicotinate](#)
 Bergamot oil, [149](#)
 Bitter orange oil, *see* [Citrus oils](#)
 Borage oil, refined, stabilised, *see* [Triglycerides](#)
 Butylenglycolum, *see* [1,3-butanediol](#)
- Cajeput, [153](#)
 Camphor oil 20%, [157](#)
 Camphor ointment, [161](#)
 Camphor spirit, [165](#)
 Caraway oil, [171](#)
 Carthami oleum raffinatum, *see* [Triglycerides](#)
 Carvi aetheroleum, *see* [Caraway oil](#)
 Caryophylli floris aetheroleum, *see* [Clove oil](#)
 Castor oil, [177](#)
 Cedar leaf oil, [191](#)
 Cedar wood oil, [195](#)
 Cera liquida, *see* [2-ethylhexyl laurate](#)
 Cetiol® HE, *see* [Polyol fatty acid ester](#)
 Cherry syrup (conserved), *see* [Syrups](#)
 Chest elixir, *see* [Wines](#)
 Chlorhexidine digluconate solution 20%, *see* [Water and aqueous solutions](#)
 Chlorhexidini digluconatis solutio 20%, *see* [Water and aqueous solutions](#)
 Cinnamon bark oil, artificial, [199](#)
 Cinnamon bark oil, Ceylon, [203](#)
 Citric oil, *see* [Citrus oils](#)
 Citronella oil, [207](#)

- Citrus oils, [211](#)
 Clary sage oil, *see* [Lavender oil / Clary sage oil](#)
 Clotrimazole 10% Cordes® RK, [219](#)
 Clove oil, [223](#)
 Coal tar spirit SR, [229](#)
 Codliver oil (type A), *see* [Triglycerides](#)
 Collodion 4%, [235](#)
 Condurango wine, *see* [Wines](#)
 Cooling cream / Cooling cream with rose oil, [239](#)
 Cooling cream DAB, *see* [Cooling cream / Cooling cream with rose oil](#)
 Cooling cream DAB 6 (stabilised, contains rose oil), *see* [Cooling cream / Cooling cream with rose oil](#)
 Coriander oil, [267](#)
 Corn oil, refined, *see* [Triglycerides](#)
 Cremor anionicus hydrophilicus SR, *see* [Hydrophilic ointments](#)
 Cremor basalis, *see* [Base cream DAC](#)
 Cremor basalis hydrophobicus, *see* [Hydrophilic ointments](#)
 Cremor nonionicus hydrophilicus SR, *see* [Hydrophilic ointments](#)
 Cypress oil, *see* [Wood oils](#)

 D-alpha tocopherol, [271](#)
 Decyl oleate, *see* [Triglycerides](#)
 Decylis oleas, *see* [Triglycerides](#)
 Decylium oleinicum, *see* [Triglycerides](#)
 Dehydol® LS 2, *see* [Macrogol-2-laurylether](#)
 Dexpanthenol, [275](#)
 Dexpanthenol concentrate 50%, [287](#)
 Dexpanthenolum, *see* [Dexpanthenol](#)
 Dexpanthenolum 50%, *see* [Dexpanthenol concentrate 50%](#)
 Dextrorotary lactic acid 20%, [293](#)
 Diethyl ether, [297](#)
 Dimethyl sulfoxide, [303](#)
 Dimethylsulfoxidum, *see* [Dimethyl sulfoxide](#)
 Dimeticon 100 / 350 / 1000, [311](#)
 Dimeticone 100, *see* [Dimeticon 100 / 350 / 1000](#)
 Dimeticone 1000, *see* [Dimeticon 100 / 350 / 1000](#)
 Dimeticone 350, *see* [Dimeticon 100 / 350 / 1000](#)
 Dimeticonum 100, *see* [Dimeticon 100 / 350 / 1000](#)
 Dimeticonum 1000, *see* [Dimeticon 100 / 350 / 1000](#)
 Dimeticonum 350, *see* [Dimeticon 100 / 350 / 1000](#)
 Dwarf pine oil, [317](#)

 Elixir e succo liquiritiae, *see* [Wines](#)

 Emulsifying hydrophobic base gel DAC, *see* [Paraffin classification](#)
 Emulsio basalis hydrophilica, *see* [Hydrophilic base emulsion \(NRF p.25\)](#)
 Ethanol 90% pure / denatured, [323](#)
 Ethanol 96% pure / denatured, [333](#)
 Ethanolum 90% rein / vergällt, *see* [Ethanol 90% pure / denatured](#)
 Ethanolum 96% rein / vergällt, *see* [Ethanol 96% pure / denatured](#)
 Ethyl acetate, [345](#)
 Ethylis acetas, *see* [Ethyl acetate](#)
 Eucalypti aetheroleum, *see* [Eucalyptus oil / Rosmary oil](#)
 Eucalyptus oil, *see* [Eucalyptus oil / Rosmary oil](#)
 Eucalyptus oil / Rosmary oil, [349](#)
 Eugenol, [359](#)
 Eugenolum, *see* [Eugenol](#)
 Evening primrose oil, refined, *see* [Triglycerides](#)
 Eye ointment, emulsifying DAC, *see* [Paraffin classification](#)
 Eye paraffin, *see* [Paraffin classification](#)

 Fennel honey (conserved), [363](#)
 Fennel oil, bitter, [367](#)
 Flaxseed oil, native, *see* [Triglycerides](#)
 Fortified wine, *see* [Wines](#)

 Gel basis Audor®, *see* [Wines](#)
 Gel Cordes®, [371](#)
 Gelatum basalis hydrophobicum, *see* [Paraffin classification](#)
 Gelatum basalis hydrophobicum emulsificans, *see* [Paraffin classification](#)
 Geranium oil, genuine, [375](#)
 Glucose syrup, [379](#)
 Glucose syrup GE 40%, *see* [Glucose syrup](#)
 Glucosum liquidum, *see* [Glucose syrup](#)
 Glycerine ointment (conserved), [383](#)
 Glycerinum 85%, *see* [Glycerol 85%](#)
 Glycerinum anhydricum >99%, *see* [Glycerol >99%](#)
 Glycerol 85%, [393](#)
 Glycerol >99%, [387](#)
 Glyceroli unguentum (konserviert), *see* [Glycerine ointment \(conserved\)](#)
 Glycerolum 85%, *see* [Glycerol 85%](#)
 Glycolic acid 70%, [413](#)
 Grape seed oil, refined, *see* [Triglycerides](#)
 Grapefruit oil, *see* [Citrus oils](#)
 Guaiazulene 25%, water-soluble, [417](#)

 Hard fat 649, [423](#)
 Honey, *see* [Honeys](#)
 Honeys, [427](#)

- Hydrated liniment SR/DAC (NRF p.40), *see* [Liniment ointments](#)
- Hydrophilic base emulsion (NRF p.25), [433](#)
- Hydrophilic cream, non-ionic SR (NRF p.26), *see* [Hydrophilic ointments](#)
- Hydrophilic cream, non-ionic without preservatives SR, *see* [Hydrophilic ointments](#)
- Hydrophilic ointment DAB, *see* [Paraffin classification](#)
- Hydrophilic ointments, [443](#)
- Hydrophobic base gel DAC, *see* [Paraffin classification](#)
- Hydroxyaethylum salicylicum, *see* [Hydroxyethyl salicylate](#)
- Hydroxyethyl salicylate, [527](#)
- Hydroxyethylis salicylas, *see* [Hydroxyethyl salicylate](#)
- Iecoris aselli oleum A, *see* [Triglycerides](#)
- int-rac-alpha-Tocopherylis acetat, *see* [All-rac-alpha-tocopheryl acetate](#)
- Ipecacuanha syrup (NRF 19.1), *see* [Syrups](#)
- Isopropyl alcohol 100%, [531](#)
- Isopropyl alcohol 70%, [557](#)
- Isopropyl myristate, [579](#)
- Isopropyl palmitate, [585](#)
- Isopropylis myristas, *see* [Isopropyl myristate](#)
- Isopropylis palmitas, *see* [Isopropyl palmitate](#)
- Jojoba oil, native, *see* [Triglycerides](#)
- Juniper oil, [591](#)
- Juniper wood oil, *see* [Wood oils](#)
- Juniperi aetheroleum, *see* [Juniper oil](#)
- Juniperi ligni etheroleum, *see* [Wood oils](#)
- Lactic acid, [595](#)
- Lanae alcoholum unguentum, *see* [Paraffin classification](#)
- Lanae alcoholum unguentum aquosum DAB/SR, *see* [Wool alcohol cream DAB/SR](#)
- Lanette® ointment (conserved), [601](#)
- Lanolin DAB, [607](#)
- Lanolinum DAB, *see* [Lanolin DAB](#)
- Laurel oil, *see* [Triglycerides](#)
- Lauromacrogol 400, *see* [Lauromacrogols](#)
- Lauromacrogols, [619](#)
- Lavender oil, *see* [Lavender oil / Clary sage oil](#)
- Lavender oil / Clary sage oil, [641](#)
- Lavendula officinalis oleum, *see* [Lavender oil / Clary sage oil](#)
- Lecithinum e sojæ, liquidus, *see* [Triglycerides](#)
- Limonis aetheroleum, *see* [Citrus oils](#)
- Lini oleum virginale, *see* [Triglycerides](#)
- Liniment ointments, [651](#)
- Linimentum aquosum N (NRF S. 39), *see* [Liniment ointments](#)
- Linimentum aquosum SR/DAC (NRF S.40), *see* [Liniment ointments](#)
- Linimentum nonionicum aquosum DAC, *see* [Liniment ointments](#)
- Liquiritiae succi elixir, *see* [Wines](#)
- Liquiritiae succus depur. Solut., *see* [Wines](#)
- Liquor aluminii acetico-tartarici, *see* [Water and aqueous solutions](#)
- Liquorice juice, purified, *see* [Wines](#)
- Lithanthracis picic spiritus SR, *see* [Coal tar spirit SR](#)
- Lotio Cordes®, [675](#)
- Macrogol 20 glycerol monostearate type II, *see* [Lauromacrogols](#)
- Macrogol 20 glyceroli monostearas, *see* [Lauromacrogols](#)
- Macrogol 4 lauryl ether, *see* [Lauromacrogols](#)
- Macrogol 8 stearate type I, [679](#)
- Macrogol ointment, [683](#)
- Macrogol-2-laurylether, [689](#)
- Macrogolglycerol hydroxystearate 40, *see* [Lauromacrogols](#)
- Macrogolglycerol ricinoleate, *see* [Polysorbate 80 / Macrogolglycerolricinoleat](#)
- Macrogolglyceroli hydroxystearas, *see* [Lauromacrogols](#)
- Macrogolglyceroli ricinoleas, *see* [Polysorbate 80 / Macrogolglycerolricinoleat](#)
- Macrogoli aether laurilicus, *see* [Lauromacrogols](#)
- Macrogoli aether laurilicus 9, *see* [Lauromacrogols](#)
- Macrogoli unguentum DAC, *see* [Macrogol ointment](#)
- Macrogols 300, [693](#)
- Macrogols 400, [697](#)
- Macrogolum 400, *see* [Macrogols 400](#)
- Mandarin oil, [705](#)
- Marshmallow syrup (conserved), *see* [Syrups](#)
- Matricaria oil, rich in (-)-alpha-bisabolol, [709](#)
- Mel, *see* [Honeys](#)
- Mel foeniculi (konserviert), *see* [Fennel honey \(conserved\)](#)
- Mel rosatum, *see* [Honeys](#)
- Melaleuca alternifolia oleum, *see* [Tea tree oil](#)
- Melissa oil, [713](#)
- Menthae arvensis aetheroleum partim mentholum depletum, *see* [Mint oil / Peppermint oil](#)
- Menthae crispae aetheroleum, *see* [Spearmint oil](#)
- Menthae piperitae aetheroleum, *see* [Mint oil / Peppermint oil](#)
- Methyl salicylate, [719](#)
- Methylis salicylas, *see* [Methyl salicylate](#)
- Methylum salicylicum, *see* [Methyl salicylate](#)

Metronidazole 25% in white petrolatum, [725](#)
 Mint oil, *see* [Mint oil / Peppermint oil](#)
 Mint oil / Peppermint oil, [729](#)
 Myrtle oil, [737](#)

Natrii laurylethersulfatis solutio 27 %, *see*
[Sodium laureth sulfate solution 27%](#)
 Natrium lacticum solutum 50%, *see* [Sodium lactate solution 50%](#)

Non-ionic hydrophil ointment, [741](#)
 Non-ionic hydrophil ointment, conserved, *see*
[Non-ionic hydrophil ointment](#)
 Non-ionic hydrous liniment SR/DAC (NRF
 P. 39), *see* [Liniment ointments](#)

Octyldodecanol, [749](#)
 Octyldodecanolum, *see* [Octyldodecanol](#)
 Oculentum emulsificans DAC, *see* [Paraffin classification](#)

Oculentum simplex DAC, *see* [Paraffin classification](#)

Oenotherae oleum raffinatum, *see*
[Triglycerides](#)

Oleic acid, [759](#)
 Oleum absinthii, *see* [Wormwood oil](#)
 Oleum amygdalarum raffinatum, *see*
[Triglycerides](#)

Oleum angelicae (e radice), *see* [Angelica root oil](#)

Oleum anisi, *see* [Anise oil / Star aniseed oil](#)
 Oleum anisi stellati, *see* [Anise oil / Star aniseed oil](#)

Oleum arachidis raffinatum, *see* [Triglycerides](#)
 Oleum aurantii dulcis, *see* [Citrus oils](#)
 Oleum aurantii pericarpium, *see* [Citrus oils](#)
 Oleum bergamottae, *see* [Bergamot oil](#)
 Oleum boraginis raffinatum, *see* [Triglycerides](#)
 Oleum Cajeputi artificiale, *see* [Cajeput](#)
 Oleum calami rectific., *see* [Acorus oil](#)
 Oleum camphoratum 20%, *see* [Camphor oil 20%](#)

Oleum cananga odorata, *see* [Ylang-ylang oil](#)
 Oleum carthami, *see* [Triglycerides](#)
 Oleum carvi, *see* [Caraway oil](#)
 Oleum caryophylli, *see* [Clove oil](#)
 Oleum cedri e fol. aether., *see* [Cedar leaf oil](#)
 Oleum cedri ligni aether., *see* [Cedar wood oil](#)
 Oleum chamomillae, *see* [Matricaria oil, rich in \(-\)-alpha-bisabolol](#)

Oleum Cinnamomi artificiale, *see* [Cinnamon bark oil, artificial](#)

Oleum cinnamomi ceylanici, *see* [Cinnamon bark oil, Ceylon](#)

Oleum citri, *see* [Citrus oils](#)
 Oleum citri paradisi artificiale, *see* [Citrus oils](#)
 Oleum citrus bergamia, *see* [Bergamot oil](#)
 Oleum coriandri, *see* [Coriander oil](#)
 Oleum cucurbita semen, *see* [Triglycerides](#)

Oleum cupressi, *see* [Wood oils](#)

Oleum cymbopogon winterianus, *see*
[Citronella oil](#)

Oleum eucalypti, *see* [Eucalyptus oil / Rosemary oil](#)

Oleum foeniculi amari, *see* [Fennel oil, bitter](#)

Oleum geranii verum, *see* [Geranium oil, genuine](#)

Oleum helianthi raffinatum, *see* [Triglycerides](#)

Oleum jecoris (Typ A), *see* [Triglycerides](#)

Oleum juglandis raffinatum, *see* [Triglycerides](#)

Oleum juniperi e baccarum, *see* [Juniper oil](#)

Oleum Juniperi e Ligno pro uso externo, *see*
[Wood oils](#)

Oleum lauri expressum, *see* [Triglycerides](#)

Oleum lavendula officinalis, *see* [Lavender oil / Clary sage oil](#)

Oleum lini, *see* [Triglycerides](#)

Oleum Mandarinum, *see* [Mandarin oil](#)

Oleum maydis raffinatum, *see* [Triglycerides](#)

Oleum melaleuca alternifolia, *see* [Tea tree oil](#)

Oleum melissae rectific., *see* [Melissa oil](#)

Oleum Menthae arvens. var. piper. japonisch,
see [Mint oil / Peppermint oil](#)

Oleum menthae arvensis, *see* [Mint oil / Peppermint oil](#)

Oleum menthae crispae, *see* [Spearment oil](#)

Oleum menthae piperitae dopp. rectific., *see*
[Mint oil / Peppermint oil](#)

Oleum myrte, *see* [Myrtle oil](#)

Oleum neutrale, *see* [Triglycerides, medium-chain](#)

Oleum oenotherae raffinatum, *see*
[Triglycerides](#)

Oleum olivarum nativum, *see* [Triglycerides](#)

Oleum omega-3 18/12, *see* [Triglycerides](#)

Oleum origani cretici, *see* [Oregano oil](#)

Oleum persea americana raffinatum, *see*
[Triglycerides](#)

Oleum petrae album rectific., *see* [Petroleum oil, white, cohobated](#)

Oleum pini pumilionis, *see* [Dwarf pine oil](#)

Oleum pini sibiricum, *see* [Spruce needle oil](#)

Oleum pini silvestris, *see* [Pine silvestris oil](#)

Oleum pruni armeniaca raffinatum, *see*
[Triglycerides](#)

Oleum rapae raffinatum, *see* [Triglycerides](#)

Oleum ricini, *see* [Castor oil](#)

Oleum rosae (verum / artificiale), *see* [Rose oil \(genuine / artificial\)](#)

Oleum rosmarini, *see* [Eucalyptus oil / Rosemary oil](#)

Oleum salviae, *see* [Sage oil](#)

Oleum salviae sclarea, *see* [Lavender oil / Clary sage oil](#)

Oleum sesami raffinatum, *see* [Triglycerides](#)

Oleum simmondsia chinensis, *see*

Triglycerides

Oleum sojæ raffinatum, *see* [Triglycerides](#)
 Oleum terebinthinae rectificat., *see* [Wood oils](#)
 Oleum thymi, *see* [Thyme oil](#)
 Oleum tritici raffinatum, *see* [Triglycerides](#)
 Oleum verbenae, *see* [Verbena oil](#),
 [nature-identical](#)
 Oleum vitis viniferae raffinatum, *see* [Triglycerides](#)
 Oleum ylang-ylang, *see* [Ylang-ylang oil](#)
 Oleyl oleate (stabilised), *see* [Triglycerides](#)
 Oleylis oleas, *see* [Triglycerides](#)
 Oleylium oleinicum, *see* [Triglycerides](#)
 Olivae oleum virginale, *see* [Triglycerides](#)
 Olive oil, native, *see* [Triglycerides](#)
 Omega-3-fish oil, *see* [Triglycerides](#)
 Oregano oil, [763](#)

Pantothenolum 50%, *see* [Dexpantenol concentrate 50%](#)
 Paraffin classification, [767](#)
 Paraffin, liquid, *see* [Paraffin classification](#)
 Paraffin, white, *see* [Paraffin classification](#)
 Paraffin, yellow, *see* [Paraffin classification](#)
 Paraffinum, *see* [Paraffin classification](#)
 Pasta Cordes®, [851](#)
 Peanut oil, refined, *see* [Triglycerides](#)
 Peppermint oil, *see* [Mint oil / Peppermint oil](#)
 Permethrin 25% solution, [857](#)
 Petroleum album rectific., *see* [Petroleum oil](#),
 [white, cohobated](#)
 Petroleum ether, [863](#)
 Petroleum oil, white, cohobated, [877](#)
 Phytomenadione, [881](#)
 Phytomenadionum, *see* [Phytomenadione](#)
 Piceae aetheroleum, *see* [Spruce needle oil](#)
 Pine silvestris oil, [885](#)
 Pine tar, [889](#)
 Pini pumilionis aetherolum, *see* [Dwarf pine oil](#)
 Pix liquida, *see* [Pine tar](#)
 Polidocanolum, *see* [Lauromacrogols](#)
 Polihexanidi solutio 20%, *see* [Polyhexanide solution 20%](#)
 Poloxyaethylenglycerinum triricinoleinicum
 35, *see* [Polysorbate 80 / Macrogolglycerolricinoleat](#)
 Polyaethylenglycol-8-stearat, *see* [Macrogol 8 stearate type I](#)
 Polyaethylenglycolum 300, *see* [Macrogols 300](#)
 Polyaethylenglycolum 400, *see* [Macrogols 400](#)
 Polyethylenglycolum laurylicum, *see* [Lauromacrogols](#)
 Polyhexanide solution 20%, [893](#)
 Polyol fatty acid ester, [899](#)
 Polyoxyaethylenglycerinum monostearinicum,
 see [Lauromacrogols](#)

Polyoxyaethylenglycerinum trioxystearinicum,
 see [Lauromacrogols](#)
 Polysorbate 20, *see* [Lauromacrogols](#)
 Polysorbate 60, *see* [Lauromacrogols](#)
 Polysorbate 80, *see* [Polysorbate 80 / Macrogolglycerolricinoleat](#)
 Polysorbate 80 / Macrogolglycerolricinoleat,
 [903](#)
 Polysorbatum 20, *see* [Lauromacrogols](#)
 Polysorbatum 60, *see* [Lauromacrogols](#)
 Polysorbatum 80, *see* [Polysorbate 80 / Macrogolglycerolricinoleat](#)
 Potash soap DAC, [909](#)
 Propylene glycol, [913](#)
 Propylenglycolum, *see* [Propylene glycol](#)
 Propylum myristinicum, *see* [Isopropyl myristate](#)
 Propylum palmitinicum, *see* [Isopropyl palmitate](#)
 Psorisol base cream Audor®, [929](#)
 Pumpkin seed oil, pressed, *see* [Triglycerides](#)

Rapae oleum raffinatum, *see* [Triglycerides](#)
 Rapeseed oil, refined, *see* [Triglycerides](#)
 Raspberry syrup (conserved), *see* [Syrups](#)
 Ribwort plantain syrup (conserved), *see* [Syrups](#)
 Ricini oleum, *see* [Castor oil](#)
 Rose honey, *see* [Honeys](#)
 Rose oil (genuine / artificial), [933](#)
 Rosemary oil, *see* [Eucalyptus oil / Rosmary oil](#)
 Rosmarini aetheroleum, *see* [Eucalyptus oil / Rosmary oil](#)

Safflower oil, refined type I, *see* [Triglycerides](#)
 Sage oil, [939](#)
 Salicylic acid 50% (paraffin), [943](#)
 Salvia officinalis oleum, *see* [Sage oil](#)
 Sapo kalinus DAC, *see* [Potash soap DAC](#)
 Sesame oil, refined, *see* [Triglycerides](#)
 Shea butter, refined, [953](#)
 Sirupus althaeae (konserviert), *see* [Syrups](#)
 Sirupus cerasi (konserviert), *see* [Syrups](#)
 Sirupus Glucosi, *see* [Glucose syrup](#)
 Sirupus ipecacuanhae, *see* [Syrups](#)
 Sirupus plantaginis lanceolatae, *see* [Syrups](#)
 Sirupus rubi idaei, *see* [Syrups](#)
 Sirupus simplex, *see* [Syrups](#)
 Skin care ointment W/L SR, *see* [Hydrophilic ointments](#)
 Skin care ointment W/L SR (without preservatives), *see* [Hydrophilic ointments](#)
 Sodium lactate solution 50%, [959](#)
 Sodium laureth sulfate solution 27%, [965](#)
 Sodium laureth sulphate solution 27%, *see* [Sodium laureth sulfate solution 27%](#)

- Soft ointment, [969](#)
 Softisan[®] 378, [981](#)
 Solutio benzalkonii chlorati 50%, *see* [Benzalkonium chloride solution 50%](#)
 Solutio Camphorae spirituosa, *see* [Camphor spirit](#)
 Solutio chlorhexidini gluconici 20%, *see* [Water and aqueous solutions](#)
 Solutio natrii laurylaethersulfurici 27 %, *see* [Sodium laureth sulfate solution 27%](#)
 Solutio polyhexanidi 20%, *see* [Polyhexanide solution 20%](#)
 Sorbitan trioleate, [991](#)
 Sorbitan- and glycerol monooleate, [987](#)
 Sorbitani trioleas, *see* [Sorbitan trioleate](#)
 Sorbitol solution 70% (non-crystallizing), *see* [Syrups](#)
 Sorbitolum liquidum non cristallisabile, *see* [Syrups](#)
 Soy lecithin, liquid, *see* [Triglycerides](#)
 Soy oil, refined, *see* [Triglycerides](#)
 Spearmint oil, [995](#)
 Spirit of ether, [999](#)
 Spirit soap, [1003](#)
 Spiritus aethereus, *see* [Spirit of ether](#)
 Spiritus camphoratus, *see* [Camphor spirit](#)
 Spiritus picis lithanthracis SR, *see* [Coal tar spirit SR](#)
 Spiritus saponatus, *see* [Spirit soap](#)
 Spruce needle oil, [1007](#)
 Squalane, [1011](#)
 Squalanum, *see* [Squalane](#)
 Star anise oil, *see* [Anise oil / Star aniseed oil](#)
 Succus liquiritiae depur. solut., *see* [Wines](#)
 Sugar syrup DAB, *see* [Syrups](#)
 Sunflower oil, refined, *see* [Triglycerides](#)
 Sweet orange oil, *see* [Citrus oils](#)
 Syrups, [1015](#)

 Tea tree oil, [1031](#)
 Thyme oil, [1035](#)
 Tocopheryl acetatis, *see* [All-rac-alpha-tocopheryl acetate](#)
 Transdermal base cream, *see* [Hydrophilic ointments](#)
 Triglycerida saturada media, *see* [Triglycerides, medium-chain](#)
 Triglycerides, [1041](#)
 Triglycerides, medium-chain, [1105](#)
 Triglycerinum diisostearanicum, *see* [Triglycerol diisostearate](#)
 Triglycerol diisostearate, [1117](#)
 Triglyceroli diisostearas, *see* [Triglycerol diisostearate](#)
 Tritici aestivi oleum raffinatum, *see* [Triglycerides](#)
 Turpentine oil, *see* [Wood oils](#)

 Unguentum adipis lanæ DAB, *see* [Lanolin DAB](#)
 Unguentum alcoholum lanæ, *see* [Paraffin classification](#)
 Unguentum alcoholum lanæ aquosum DAB/SR, *see* [Wool alcohol cream DAB/SR](#)
 Unguentum alcoholum lanæ SR, *see* [Paraffin classification](#)
 Unguentum basale, *see* [Base cream DAC](#)
 Unguentum basale hydrophobicum, *see* [Hydrophilic ointments](#)
 Unguentum camphoratum, *see* [Camphor ointment](#)
 Unguentum Cordes[®], *see* [Paraffin classification](#)
 Unguentum emulsificans aquosum cum propylenglycolum, *see* [Hydrophilic ointments](#)
 Unguentum emulsificans aquosum DAB, *see* [Hydrophilic ointments](#)
 Unguentum emulsificans aquosum N sine conservante SR, *see* [Hydrophilic ointments](#)
 Unguentum emulsificans aquosum N SR (NRF S. 26), *see* [Hydrophilic ointments](#)
 Unguentum emulsificans aquosum sine conservante, *see* [Hydrophilic ointments](#)
 Unguentum emulsificans aquosum SR (NRF S. 27), *see* [Hydrophilic ointments](#)
 Unguentum emulsificans DAB, *see* [Paraffin classification](#)
 Unguentum emulsificans nonionicum aquosum DAB, *see* [Non-ionic hydrophil ointment](#)
 Unguentum glycerini (konserviert), *see* [Glycerine ointment \(conserved\)](#)
 Unguentum Lanette[®], *see* [Lanette[®] ointment \(conserved\)](#)
 Unguentum leniens, *see* [Cooling cream / Cooling cream with rose oil](#)
 Unguentum leniens DAB 6, *see* [Cooling cream / Cooling cream with rose oil](#)
 Unguentum molle, *see* [Soft ointment](#)
 Unguentum ophthalmicum emulsificans DAC, *see* [Paraffin classification](#)
 Unguentum polyäthylenglycoli, *see* [Macrogol ointment](#)

 Vaseline, white, *see* [Paraffin classification](#)
 Vaseline, yellow, *see* [Paraffin classification](#)
 Vaselinum album, *see* [Paraffin classification](#)
 Vaselinum album pro ophthalmicum, *see* [Paraffin classification](#)
 Vaselinum flavum, *see* [Paraffin classification](#)

Verbena oil, nature-identical, [1121](#)
Vinum condurangum, *see* [Wines](#)
Viscous base solution (NRF p. 20), [1125](#)
Walnut oil, refined, *see* [Triglycerides](#)
Water, *see* [Water and aqueous solutions](#)
Water and aqueous solutions, [1137](#)
Wheatgerm oil, refined, *see* [Triglycerides](#)
Wines, [1161](#)
Wolff® base half-fat cream, *see* [Hydrophilic ointments](#)
Wood oils, [1167](#)

Wool alcohol cream DAB/SR, [1173](#)
Wool alcohol ointment DAB, *see* [Paraffin classification](#)
Wool alcohol ointment SR with white paraffin, *see* [Paraffin classification](#)
Wool alcohol ointment SR/DAC, *see* [Paraffin classification](#)
Wool fat, [1201](#)
Wormwood oil, [1215](#)

Ylang-ylang oil, [1219](#)