

Validation documentation
Narcotic substances, solid

HiperScan GmbH

December 6, 2019

Contents

Introduction	5
Context of this document	5
Criteria for the inclusion of substances	6
Validation concept	7
Model creation procedure and validation runs	7
Particularities of individual substance classes	9
Significance of testing with <i>Apo-Ident</i>	10
Conclusion	11
Explanation of terminology	12
Summary	17
Validation samples	17
Validation results	18
Validation reports	19
Amphetamine sulfate	19
Barbital sodium	25
Cocaine hydrochloride	29
Codeine phosphate hemihydrate	33
Dexamphetamine sulfate	37
Dihydrocodeine hydrogen tartrate	41
Fentanyl citrate	45
Hydromorphone hydrochloride	49
Methadone hydrochloride	53
Midazolam	61
Midazolam hydrochloride	65
Morphine hydrochloride	69
Morphine sulfate	75
Oxycodone hydrochloride	79
Phenobarbital	83
Appendix	87
Additional calibration samples (<i>Type A</i>)	87
Additional validation samples (<i>Type B</i>)	125
Additional validation samples (<i>Type C</i>)	211
Requirements of validation	729
Conformity of <i>Apo-Ident</i> with the <i>European pharmacopeia</i>	730
Addition to the two-stage models	731
References	733
Index	735

Introduction

The unequivocal proof of the identity of pharmaceutical raw materials based on a monograph or traditional alternative methods is work-intensive, time-intensive and economically often no longer makes sense. Near-infrared spectroscopy (NIR) provides a new means here. It enables relatively easy, fast and nonetheless safe identity testing via the preparation and evaluation of spectra.

The analysis system *Apo-Ident* was developed specially for pharmacy use. Pharmacists have the duty to check the identity of all raw materials for extemporaneous products in their pharmacies. This normally takes place based on the monographs for the respective substances in the European Pharmacopoeia. But also NIR spectroscopy is described in the European Pharmacopoeia as an identification method which, as opposed to the methods incorporated in the respective monographs, is approved for testing [1] (quotation translated)

with the prerequisite that the same results (“namely the determination of the identity” [2]), (quotation translated) are achieved as with the described methods and instruments.

The *Apo-Ident* analysis system serves to identify raw materials for prescriptions in the defined manner according to *ApBetrO* [Pharmacies Rules and Regulations] §§ 6 and 11 at pharmacies (NIR spectroscopy as an alternative testing method). *Apo-Ident* consists of three components:

- An *NIR spectrometer*, which records the spectra of non-preprocessed raw materials in a measuring glass in diffuse reflection or transfection.
- The *QuickStep* spectroscopy software controls the instrument and records the spectra and user inputs via a pharmacy-specific software plug-in. It also generates the test protocol for documentation of the testing and storage of the printout to be signed in the pharmacy.
- The software module *IdentModule* incorporates *reference databases*. The spectra from the *QuickStep* software are presented to it for evaluation.

NIR spectroscopy is a very powerful analytical method. It is also able to establish the identity of several chemical compounds and mixes in as far as an appropriate database (technically correct: a [chemometric model](#)) was created. Identity testing with *Apo-Ident* is a very safe, very fast and easy to operate analytical method for testing a large number of raw materials.

Context of this document

The suitability of the instrument, method and database is proven as follows:

- *NIR spectroscopy as a method for identity testing*: The *Ph. Eur.* [3] describes NIR spectroscopy in *Section 2.2.40* as an analytical method which is also suitable for the identification of raw materials. Therefore, validation of the method as such is not necessary.
- *Performance of the instrument*: The *Ph. Eur.* [3] furthermore describes the apparatus and the testing of its performance in *Section 2.2.40*. The document *Erfüllung von 2.2.40 Ph. Eur. durch Apo-Ident* [4] compares the implementation by *Apo-Ident* with this monograph in order to prove that *Apo-Ident* meets the specifications of the Pharmacopoeia. Each individual instrument delivered to a pharmacy is qualified in accordance with the tests described in “*Control of Instrument Performance*”. In this test, the unit consisting of analysis instrument hardware and the *QuickStep* spectroscopy software is assessed. The result is documented in a test protocol which is kept at the pharmacy.
- *Validation of the database* is documented separately for each substance class. The report at hand documents the substance class *Narcotic substances, solid*.

The *Arbeitsgemeinschaft der Pharmazieräte Deutschlands (APD)* [Working Group of German Pharmacy Inspectors] has clarified the following in its resolution dated October 16, 2013 ([5], quotation translated):

NIR is a testing method incorporated in the Pharmacopoeia. The testing quality depends on the quality of the database stored. The APD views the use of NIR instruments in case of ensured validation of the databases used in conjunction with it as one of several options for identity testing.

The APD defined more precisely ([6], quotation translated) on October 1, 2014:

The use of near-infrared is a recognised testing method according to Ph. Eur. 8. For the use of NIR instruments in pharmacies for testing the identity of raw materials, sufficient and verifiable validation of the instrument used is required. The quality of the database stored by the instrument manufacturer is decisive for quality. Batch-specific differences with the same original substances must be taken into account if present.

So NIR is basically suitable. The validity of the reference database is proven with the existing validation documentation.

Criteria for the inclusion of substances

This validation documentation describes the results of the validation of the reference database for the substance class *Narcotic substances, solid*. Validation documentation is created for each published version of the reference database for all substance classes incorporated.

The reference database is incorporated in the software module *IdentModule*. During identification testing with *Apo-Ident*, spectra which are used for evaluation purposes are presented to it by the *QuickStep* software. In the same manner, the *IdentModule* is presented all validation spectra successively during the validation runs for evaluation purposes. The *IdentModule* responds respectively (without taking the initial assumption into account) with the identified substance or rejects it as unknown. The correctness of this response is checked for each possible initial assumption and counted.

The results are summarised for each substance and reproduced in this document. The core statement of this validation report is that the following criteria must be fulfilled for each database entry, so that *Apo-Ident* offers verification of identity by means of NIR for the relevant substance/substance group:

- The database is exclusively generated from spectra which have been recorded by *HiperScan GmbH* on traceable samples in pharmaceutical quality.
 - The samples are procured via typical pharmacy sources (*DAC III.2.: Bezugsquellennachweis für Rezepturbestandteile [reference source for prescription components]* [7]).
 - A valid manufacturer's certificate exists (content, purity and identity of the batch).
 - The identity was confirmed by a certified test laboratory or *HiperScan GmbH*.
- Each version of the reference database (every update) is validated in-full.
 - Calibration spectra (*Type A*), other spectra recorded under the control of *HiperScan GmbH* (*Type B*), and spectra from the field (*Type C*) are presented to the *IdentModule* for evaluation in three separately evaluated validation runs.
 - Here, no single *false positive* result may arise.
 - Here, the various substance classes are also tested for reciprocal rejection, where this is objectively justified (see *Summary* section).
- In the validation with spectra recorded under the control of *HiperScan GmbH*, spectra of at least one independent sample must be considered, i.e. spectra from a batch whose spectra have not been used for the generation of the database. In addition, the set of *Type A* and *Type B* spectra must originate from at least three different batches.
- Spectra of additional substances may be used for the generation of the database even though they will not be offered for identification by this database. The purpose is the reliable distinction from these substances.

- Any positive result of *Apo-Ident* confirms the identity of the substance/substance group and distinguishes it from all other substances of the database. In the case of substance groups, the result is ambiguous: Distinction from all substances not belonging to the group is proven. The substance is identified as a member of this group. However, within the substance group, it is not possible to reliably classify which substance has been tested.
- The criteria for clear identifiability are a [specificity](#) of 100 % ([true negative rate](#)) and a minimum distance in the distance matrix. See 2. d) under [Model creation procedure and validation runs](#).

Validation concept

[Chemometrics](#) is a statistical technique for the extraction of relevant chemical information from spectra. In mathematics, this method is described as *multivariate data analysis*. Chemometrics proceeds here as follows:

1. Collection of spectra for the *calibration sample*. The results (identities) of the calibration sample must be known. The calibration samples must be representative for the samples which are to be evaluated later. Therefore, they must take the various possible (physical) compositions into account. (Therefore, sourcing calibration samples for NIR from the specialist trade is superior to the use of CRS reference substances.)
2. The first mathematical step is *calibration*. Here, the [chemometric model](#) is calculated from the *calibration sample spectra* ([reference spectra](#)) and limits as well as some parameters are stipulated. The chemometric model is used later to calculate the analysis result (*prediction*).
3. Collection of further spectra for the *validation sample* which should be independent of the *calibration sample*. The results (identities) of the *validation sample* must also be known. The textbook suggests a random sample with a normal scope of 25 % to 50 % of the *calibration sample* [8].
4. The second technical data step is *validation*. Here, the [chemometric model](#) created is evaluated based on the spectra of the *validation samples*. As validation parameters for the identification, the *Ph. Eur. Section 2.2.40* [3] specifies the [specificity](#) and [robustness](#).

The validation step according to the textbook has the target of estimating the performance capability of the model created based on a random sample. In order to achieve the best possible precision, attention is paid to the calibration sample. In the field of pharmaceuticals, the safety of the method has priority. In order to be able to [validate](#) the model within the regulatory scope, the validation step must include probative force. For this purpose, the validation sample must be *representative and complete* in order to enable the testing of all cases.

A *sufficient number of batches* must be secured for validation because validation finally proves whether the number of batches in calibration suffices.

Each substance is validated individually. The validation results are documented per substance in this document. Moreover, the documents show how many and which batches have been used for creating the model or model validation.

At least one certificate is taken in for each substance from an accredited test laboratory for the independent testing of identity of the sample. The identification number of the corresponding test certificate is listed in the report, enabling traceability of a substance tested according to the monographs in the Pharmacopoeia.

Model creation procedure and validation runs

The safety of the [chemometric models](#) is guaranteed by several measures during model creation, of which the validation step is the final one. Normally, the procedure is as follows. It is in particular valid for the active pharmaceutical ingredients (APIs) *solid API excipients*, *liquid/semi-solid API excipients (with a test certificate)*, *narcotics - solid medicinal substances* and *drugs*. If, for individual substance classes, variations are required, they are depicted in the section [Particularities of individual substance classes](#).

1. Collecting the reference spectra (calibration sample)

- a) Procurement of the samples from the same sources from which pharmacies source their raw material for compounding (Caelo, Fagron, Euro-OTC, . . . , see also *DAC III.2. Bezugsquellen-nachweis für Rezepturbestandteile* [Sources of supply for compounding] [7]).
- b) Testing the suitability according to *ApBetrO* [Pharmacies Rules and Regulations] §§ 6, 11, that is to say the availability of a valid manufacturer certificates via identity, purity and contents of the batch.
- c) Recording standard 40 spectra of the sample in different positions, as a standard on four instruments. Here, handling and presentation of the samples as later in the pharmacy.
- d) Visual checking for anomalies in the spectra. In case of indications of measurement errors, measurement must be repeated. If a signature is missing in the spectrum, the substance may be excluded from the start as not promising (the spectra are nonetheless entered in the database validation as independent *Type B* spectra).
- e) Testing identity. For each substance, a certificate of correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the respective following substance page of this validation documentation the *Mahalanobis distance* to this reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance. Such samples underline the statistical spread of the original reference substance, but cannot add any new characteristics of the substance.

HiperScan GmbH cooperates with some suppliers as follows: the raw materials supplier takes a sufficiently large sample in his incoming goods area so that part of it can be used for recording the NIR spectra. The remainder of the sample goes to analytics for market approval. The manufacturer's batch certificate emerges from these identity, contents and purity tests, which consequently also prove the correct identity of the NIR reference sample. Therefore, the NIR spectra are suitable for structuring the database (*Type A*) and can also optionally be used for validation (*Type B*). The samples which this applies to are marked in the validation report with a footnote.

- f) If the identity of the new sample is proven, it is declared as a reference sample and the spectra are approved for structuring the database.

2. Generating the chemometric models (calibration)

- a) Determination of the transformation matrix from the reference spectra using variance maximisation [8, 9]. (*All* reference spectra are always included, even if only a few spectra are added for an update.) All reference spectra receive the same data pre-treatment, which is also later applied to all measurement spectra in the field (in the pharmacies).
- b) Checking that the number of principal components used is still sufficient.
- c) Calculating the limits for each substance from the spread of the reference spectra. The calculation regulation is identical for each substance in a substance class.
- d) Checking the distances between the limits of the separable substances: the distance matrix incorporates the *Mahalanobis distances* from each substance to all the others. If a distance is less than 10, both two substances must be declared as inseparable or *both two* must be removed from the database. (The spectra of removed substances remain as independent *Type B* spectra in the validation.)
- e) Testing the model based on the reference spectra. If *false positive* results arise, you must proceed as in the case of *Mahalanobis distances* which are too small.
- f) If chemometric models are available for all substance classes which meet both criteria (distance matrix and no *false positives*), they are joined together with the evaluation algorithms as an *IdentModule* and encrypted. This unit can no longer be changed. Its overall function is tested by validation.

3. Set of validation spectra (validation samples)

The following is provided for validation:

- a) *Type A*: The reference spectra = calibration spectra from which the database was generated. These also include spectra from substances which the *chemometric model* should not identify, but were also recorded during generation in order to increase selectivity. (As a result, the model “learns” to differentiate from other substances which are actually unknown to it.)
- b) *Type B*: Spectra recorded under the control of *HiperScan GmbH* and not used for the generation of the database. These also include reference spectra of other substance classes, and spectra that are not used as reference spectra. Samples are considered to be independent, if they originate from a batch, of which no spectra have been used for the generation of the database. (Up to *IdentModule 2018-01*, samples were still considered independent if the sampling was done independently, i.e. if they originate from another sales container.)
- c) *Type C*: Spectra from the field, which have not been recorded under the control of *HiperScan GmbH* and have not been used for the generation of the database. The spectra include not only substances of the substance class to be tested, but also substances from other classes.

All manufacturers’ batches from which spectra are used for the validation are listed by substance in this document: for substances included in the substance class *Narcotic substances, solid* in the respective validation reports; otherwise in attachments *A*, *B* and *C*.

Furthermore remains valid: validation spectra may only be removed if a error in the spectrum can be proven. Here, the spectra are not deleted, but instead placed on a *blacklist* incorporating the reason, date and initials in the commentary.

The section *Particularities of individual substance classes* treats the other substance classes from which *Type B* and *Type C* spectra are cited for validation purposes.

4. Validation runs and approval

- a) Validation spectra are transferred holistically to the *IdentModule* for evaluation in the same way as the spectroscopy software *QuickStep* transfers measured spectra.
- b) Following the provision of each spectrum, the *IdentModule* responds as to whether it has recognised a substance and which substance was recognised.
- c) The correctness of this response is checked for each possible initial assumption (each measurable substance with the substance class) and counted according to *true negative*, *false negative*, *true positive* and *false positive*. These figures are provided for each substance and additionally in the section *Summary*, separated according to types *A*, *B* and *C*.
- d) No *false positive* results whatsoever are permissible.
- e) If the criterion is also met for all substance classes, the *IdentModule* is approved.

Particularities of individual substance classes

Basically, *HiperScan GmbH* procures and tests the manufacturer’s certificate for the batch, commissions external testing of the identity of the sample or carries it out independently and stores the certificates. As described, this process is established for the Pharmacopoeia substances, that is to say for substance classes **APIs & excipients, solid**, **APIs & excipients, liquid/semi-solid (with analysis certificate)**, **Narcotic substances, liquid/semi-solid** and **Drugs**. Therefore, *HiperScan GmbH* is able to furnish proof of the identity of the reference samples. In case of manufacturer-specific substance classes and others, individual steps are organised differently in-part:

The substance class **APIs & excipients, liquid/semi-solid (other)** (often described as cosmetics) incorporates substances for which no specification of the requirements of the pharmaceutical quality is determined, neither in a Pharmacopoeia monograph, a DAC/NRF monograph nor via a manufacturer’s specification. Consequently, neither the identity nor contents can be tested independently. No certificates whatsoever exist for the reference samples. So here, merely the matching of the sample with former samples of this product is established and confusion with the other substances

is ruled-out. (If the manufacturer of such a substance prepares a specification, determines testing methods and provides manufacturer's certificates in accordance with *ApBetrO* [Pharmacies Rules and Regulations] §§ 6, 11, *HiperScan GmbH* can assign the substance to the substance class *APIs & excipients, liquid/semi-solid (with analysis certificate)* again in the future).

Substance class **HCK – nutritional supplements (Hepart)** contains the HCK micro-nutrients from the Swiss company *Hepart AG*. *HiperScan GmbH* receives the reference samples directly from the manufacturer. For each reference sample, *HiperScan GmbH* also receives manufacturer's certificates and keeps these. New checking of the identity of the reference sample is not carried out by *HiperScan GmbH*. The identity of the reference samples is therefore documented by *Hepart AG*. The spectra of all batches provided by *Hepart AG* are recorded by *HiperScan GmbH* and entered in the database.

All the manufacturer's batches are used for the generation and validation of the substance class *HCK*. The expected variation is also represented in the generation and validation if there are less than three batches.

Also, for the substance class **PhytoComm** (TCM-Granulated herbal extracts of the manufacturer *PhytoComm*) spectra for all useable batches are recorded by *HiperScan GmbH* and entered in the database. The supplier organises the respective tests themselves and keeps the test certificates.

A new evaluation option was created for the class *PhytoComm* with the update 2016-01. As the risks are considerably fewer than those from chemical agents, the pharmacist can specify a reasonable criterion for the *specificity* in accordance with internal risk estimation. The database for this is created without taking safety distances into account and no criterion is determined in advance for the *specificity*. Instead, the *specificity* for testing the identity with this concrete substance is calculated in the validation for each substance and provided with the measurement result. The pharmacist then judges himself whether this safety is reasonable with regard to the risk of the substance.

Additionally, a statistical forecast is provided for the *specificity* which is determined according to the *Rule of Three* [10, 11]. For this forecast, it is assumed that there would have been three wrong results more and is provided with a lower limit for *specificity*. This value has a special meaning if a *specificity* of 100 % is achieved for a substance during validation. In this case, the lower limit allows conclusions regarding the scale of existing safety for which with an endless number of validation spectra a value of less than 100 % is to be assumed.

If, for example 14000 spectra not belonging to the substance are presented and no *false positive* classification is made, a hypothetical number of three *false positive* results is assumed (*Rule of Three* [10, 11]) and the *specificity* is defined with 100.0000 % (> 99.9786 %). Here, it applies that the higher the number of validation spectra which form the statistical basis, the better the *specificity* calculated via the lower *specificity* limit will be approximated.

The positive result of the identity test using *Apo-Ident* establishes that the sample spectrum is in accordance with a batch of the specified granulate from the supplier *PhytoComm*, whereby all useable batches from the supplier are known.

The *PhytoComm* class can only confirm the identity of batches that have been used for the generation of the database. As a consequence there cannot be any validation spectra of other batches. Therefore, the criterion reads that two samples (from different sales containers) from each batch must exist, one for the structure of the database (*Type A*) and one for the validation (*Type B*).

Significance of testing with *Apo-Ident*

The analysis result is determined using sophisticated statistical methods according to state-of-the-art science and technology. Chemical and pharmaceutical knowledge is applied for the selection of the samples from which the calibration spectra and validation spectra are recorded. Otherwise it does not influence the further steps of model creation.

Verbally, the statement of the analysis result can be expressed as follows. Here "*the spectra match*" means that the criteria *Mahalanobis distance*, *outlier analysis* and *correlation* are met as shown in *Erfüllung von 2.2.40 Ph. Eur. durch Apo-Ident* [4]. "The spectra do not match", on the other hand, means that at least the criterion *Mahalanobis distance* is not met.

The positive analysis result "*was identified as ...*" is very meaningful because both the quantity of substances to be taken into account and the number of underlying samples is very comprehensive.

1. The spectrum of the sample measured matches spectra of the defined substance.

2. The spectrum of the sample measured does not match any spectrum of any other substance in this substance class. Therefore, all other substances can be clearly ruled-out.
3. As the spectra from other substance classes were used for validation, it is proven that no spectrum of one of these other substances matches the defined substance. (All substance classes with which a spectrum comparison is possible and makes sense are used for validation. This is documented for each substance class in the section [Summary](#).)
4. If the defined substance belongs to a [substance group](#) which in itself is not clearly separable with *Apo-Ident*, matching with the spectra of one or several substances in this group is confirmed. Which of these substances it actually is cannot be determined clearly. All other substances are excluded analogous to [2](#) and [3](#).

On the other hand, a negative analysis result “*was not identified as ...*” means:

1. The substance offered could not be recognised based on the spectrum of this sample.
2. The identity of this sample is not confirmed.
3. Testing must be repeated in accordance with the specifications of the Pharmacopoeia.

Conclusion

NIR spectroscopy is a testing method incorporated in the Pharmacopoeia. In case of successful database validation, it is a possible method for identity testing [\[5\]](#). *Apo-Ident* meets the criteria of the *European Pharmacopoeia* as a near-infrared spectrometer and proves the validity of the reference database with the existing validation documentation. This means that *Apo-Ident* can be used as an alternative testing method for testing raw materials at pharmacies.

Explanation of terminology

The following section serves to explain or define specialist terminology which is required in order to understand this document. If necessary, definitions for the analysis system *Apo-Ident* are defined more precisely.

The term database is used in this document exactly as in the *Ph. Eur. Section 2.2.40* [3] synonymous with [chemometric model](#). In order to differentiate the databases which are relatively independent of each other, *HiperScan GmbH* frequently also uses the term [substance class](#) (primarily in the plural). On the other hand, the spectra used to structure the database are termed spectrum collection and not database.

Substance classes are units of the organisational structure of the *IdentModule*. The substance classes are substance [databases](#) which are also broadly independently subscribable. On the one hand, the liquid and semi-solid substances are separated from the solid powders because they are measured against different references and therefore the spectra cannot be compared. On the other hand, for example the Pharmacopoeia substances are kept separated from the manufacturer-specific database *PhytoComm* for TCM (traditional Chinese medicine) raw materials.

The individual substance classes need only be limited against each other in-part. Often, no risk of confusion exists because they can only be procured from different sources. On the other hand, in several cases we handle substances which need not be distinguished. For example, em Huang Qi granulate from the company *PhytoComm* neither needs to be delimited from *Huang Qi* granulate from the company *HerbaSinica* nor is matching required. Respectively one single [chemometric model](#) is behind a substance class. (Even if several reciprocally secured chemometric models would be permissible.) The terms *substance class*, *chemometric model* and *databases* are mostly used here as synonyms.

A substance group respectively summarises all the substances within a [substance class](#) which cannot safely be distinguished from one another based on their NIR spectra. However, all the other substances in the database can be excluded.

The formation of subgroups is mentioned in the *Ph. Eur. Section 2.2.40* [3]. In this manner, technical restrictions in case of extensive databases can be avoided and it is possible to prepare individual subgroups with different spectrum pre-treatment. Validation of the subgroups against each other is required. *HiperScan GmbH* has solved these technical restrictions and doesn't use any subgroups within a substance class any longer.

Principal component analysis (PCA) [8, 9] is a multivariate statistics process or multivariate data analysis. It serves to structure, simplify and illustrate comprehensive data records by describing a large number of statistical variables by describing a lower number of linear combinations (the *principal components*) which are as significant as possible. In the *Apo-Ident IdentModule*, *PCA* is used to evaluate the recorded spectrum data (corresponding with *Ph. Eur. 2.2.40* [3]).

The term validation is defined in both relevant contexts here with different (even if related) meanings.

Within the sense of the expert discipline of *chemometrics*, validation is a process step when creating a [chemometric model](#): after a transformation matrix, limits and various parameters have been calculated or determined from a set of reference spectra during the course of the calibration step [8, 9], the validation step determines the performance capability of the model (selectivity, precision, ...) based on the validation spectra. Normally, random sampling is planned here. In order for the validation to gain strength of proof, the validation spectrum set must be selected with an appropriately wide scope (*representative* and *complete*). The terms *validation run* and *validation step* always actually mean the process step in this sense.

In the regulatory sense (of pharmaceutical production), validation is the documented proof that a process or system meets the previously specified requirements reproducibly when applied practi-

cally. In this sense, the *Apo-Ident* databases only become validated databases with the validation documentation, which this document is part of.

The *European Pharmacopeia* uses the term validation in *Section 2.2.40* within the sense of the specialist discipline of *chemometrics* [3].

The robustness of a process is the property of only being influenced by environmental fluctuations (e.g. temperature or humidity) a little. A method is robust if the environmental conditions do not or hardly falsify the final result.

The specificity of a classification (of a [chemometric model](#)) is the [true negative rate](#).

The recognition rate (also sensitivity) is the [true positive rate](#). It defines in how many percent of cases a correctly set up substance is actually confirmed.

The true negative rate describes the share of spectra correctly classified as non-identity during validation. This is equivalent to correct classification. It means that a substance *A* within identity checking as substance *B* is judged as “*not identified*”. The *true negative rate* is equivalent to the conditional frequency

$$h(\text{rejected}|\text{genuinely no identity}) = \frac{r_n}{r_n + f_p}$$

with r_n as the total number of *true negative* classifications and f_p as the total number of *false positive* classifications. For successful validation of an *IdentModule*, all spectra presented belonging to this category must be classified as *not in accordance*.

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of this number. The weight of each spectrum of a substance/substance group *i* therefore results as

$$w_i = \frac{1}{n_i}$$

with n_i number of spectra in this substance/substance group. This weighting ensures that the overall result cannot be enhanced by adding especially large numbers of spectra from easily separable substances.

The true positive rate describes the share of spectra correctly classified as identity during validation. This is equivalent to correct classification. It means that a substance *A* within identity checking as substance *A* is judged as “*identified*”. The *true positive rate* is equivalent to the conditional frequency

$$h(\text{identified}|\text{genuine identity}) = \frac{r_p}{r_p + f_n}$$

with r_p as the total number of *true positive* classifications and f_n as the total number of *false negative* classifications. The *true positive rate* is a measure for the recognition rate of the validated *Apo-Ident IdentModule*.

In order to ensure that each substance is received with the same weight, the spectra are weighted as described for the [true negative rate](#).

The true negative result describes a spectrum correctly classified as non-identity during validation. It is equivalent to correct classification. It means that a substance *A* within identity checking as substance *B* is judged as “*not identified*”.

The false positive result describes a spectrum falsely classified as non-identity during validation. This is the most critical type of possible false classification. It means that a substance *A* within identity checking as substance *B* is judged as “*identified*”. For successful validation of an *IdentModule*, a number of false positive events of zero are demanded for all spectra entering the validation. The exception to this restriction is the class of TCM granulates from the company *PhytoComm* as described under [Particularities of individual substance classes](#).

The true positive result describes a spectrum correctly classified as identity during validation. It is equivalent to correct classification. It means that a substance *A* within identity checking as substance *A* is judged as “*identified*”.

The false negative result describes a spectrum falsely classified as non-identity during validation. It is equivalent to false classification. It means that a substance *A* within identity checking as substance *A* is judged as “*not identified*”.

The ‘Rule of Three’ says that with a probability of 95 % the next random sample of the same size no more than three false results are to be expected if no false result existed in the existing random sample [10, 11].

The *specificity* and *recognition rate* are determined both globally and from the validation runs for all substances. The information is supplemented with the hypothetical value if there had been three false results more. The percent information is provided in parentheses with the “greater than” symbol ‘>’, e.g. *specificity* 100.000 % (>99.983 %) if 17 567 false spectra have been presented without one single *false positive* result. The larger the statistical basis, the lower the influence of the hypothetical false results.

The Mahalanobis distance is a distance measure between two points in *n*-dimensional vector space. Here, the respective direction component of the distance to *standard deviation* [12] of an *n*-dimensional distribution is standardised. In case of the *principal component analysis* [8, 9] this standardisation relates to the distribution of the respective calibration data set for a classification (substance/substance group) in the *principal component space* [8]. The *Mahalanobis distance* of a point (mapping of a spectrum) $\vec{y}$ in the *n*-dimensional principal component space to the expected value of an *n*-dimensional distribution $\mathbf{X}$ then results as

$$d(\mathbf{X}, \vec{y}) = \sqrt{(\bar{\mathbf{X}} - \vec{y})^T \mathbf{S}^{-1} (\bar{\mathbf{X}} - \vec{y})} \quad \text{with} \quad \mathbf{X} \in \mathbb{R}^{m \times n}, \vec{y} \in \mathbb{R}^m$$

[13]. Here, *m* is equivalent to the number of principal components used (dimension of the principal component space) and *n* the number of measurements existing in the calibration data set (spectra). $\bar{\mathbf{X}}$ is the expected value of the resulting distribution for the calibration data set (the average value of *n* measurements received). $\mathbf{S}^{-1}$ is the inverse covariance matrix [12] for distribution $\mathbf{X}$.

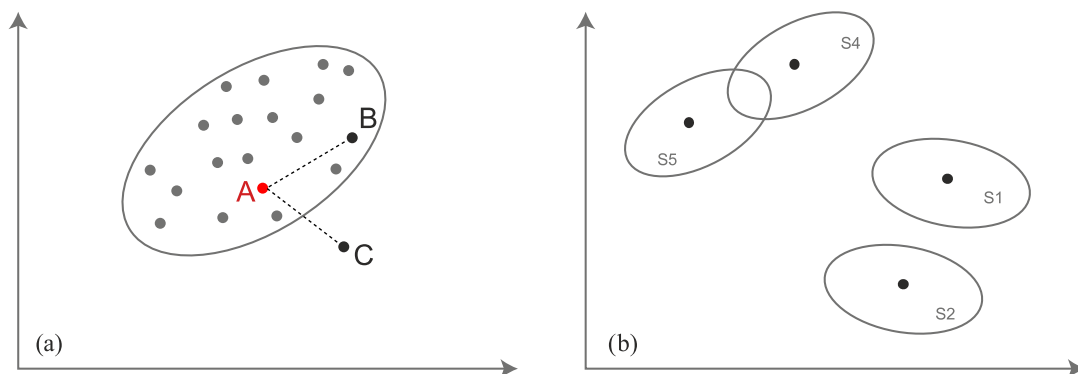


Figure 1: (a) The *Mahalanobis distance* from *A* to *B* is less than from *A* to *C*. However, the Euclidean distances are identical. (b) The *Mahalanobis distance* between the two measurement series *S4* and *S5* is smaller than between *S1* and *S2*. However, the Euclidean distances are identical.

The *Mahalanobis distance* offers advantages compared to the *Euclidean distance*: For the calculation of the distance it takes the statistical properties of a data point distribution (measurement series), i.e. average value, variance and covariance of the data points [14] into account. The *Mahalanobis distance* is applied while creating the reference database for evaluating the spectra from different samples of a substance.

A chemometric model is a classifier based on statistical methods [8, 9]. Through the respective algorithm used (e.g. *Principal Component Analysis*, *Cluster Analysis*), a maximum of chemical information is extracted from measurement data. Here, systematical or physical disturbances are eliminated using appropriate data pre-processing [15, 16].

At several places in this document, in order to simplify understanding, the term **database** is used instead of *chemometric model* – in the same manner as in the *Ph. Eur. Section 2.2.40* [3].

A sample (with its own sample ID) refers to substance in a sales container. Repeated sampling from the same sales container is listed under the same sample ID. (The suffix “SI” is not part of the sample ID.) Several samples may originate from the same batch. Samples are called “independent”, if they originate from a batch, of which no spectra have been used for the generation of the database. (Up to *IdentModule 2018-01*, samples from different sales packages were considered to be independent.) The information above the list of validation spectra now includes also the number of batches that deliver independent samples for the validation (for both *Type B* and *Type C*).

In case a supplier takes a sample for testing from its incoming goods and splits it to multiple laboratory containers, the substance in all laboratory containers will still be ascribed to the same sample. *HiperScan GmbH* only uses one of the subsamples.

Reference samples are used to structure the database. The *reference spectra* originate from these samples. In chemometric technical jargon you would normally say: For *calibration*, a **chemometric model** is generated from the *calibration spectra* recorded from the *calibration samples*, whose quality is subsequently assessed in **validation**.

Reference samples are procured via typical pharmacy sources. Their identity is tested. The *reference spectra* are recorded by *HiperScan GmbH*. The documentation also includes the manufacturer’s name and batch number.

Reference samples are clearly identified by a sample ID. Samples without sample ID may not be used as *reference samples*.

Summary

A total of 446 856 spectra from 34 779 different batches for a total of 1223 substances were used to validate the substance class *Narcotic substances, solid*.

Validation samples

The validation samples can be categorised as follows:

Type A Calibration spectra. These are the spectra used to generate the chemometric model. They were recorded by *HiperScan GmbH*. Detailed information regarding the batches or samples can be found in the following validation reports under *calibration samples* and under *Type A*. Further information is listed in [Appendix A](#).

Substance class	Substances	Batches	Spectra
Narcotic substances, solid	15	75	6951
APIs & excipients, solid	349	1789	99 518
Specific classifier 26188	10	20	980
Drugs	15	18	800

From category *A* a total of 108 249 spectra from 1902 batches for a total of 389 substances were taken into account for validation.

Type B Spectra from independent samples which are not included in database generation. These spectra were recorded by *HiperScan GmbH*. Detailed information regarding the batches or samples can be found in the following validation reports in the section *Type B* or in [Appendix B](#).

Substance class	Substances	Batches	Spectra
Narcotic substances, solid	17	58	4089
APIs & excipients, solid	505	1700	83 707
Drugs	263	490	22 603
Specific classifier 26188	10	14	589
TCM - Granulated herbal extracts (PhytoComm)	163	836	90 893
HCK - nutritional supplements (Hepart)	120	812	41 085

From category *B* a total of 242 966 spectra from 3910 batches for a total of 1078 substances were taken into account for validation.

Type C Spectra from independent samples which are not included in database generation. *Apo-Ident* customers carried out the measurements. Detailed information regarding the batches or samples can be found in the following validation reports in the section *Type C* or in [Appendix C](#).

Substance class	Substances	Batches	Spectra
Narcotic substances, solid	15	449	1003
APIs & excipients, solid	397	23 086	78 141
Drugs	284	5792	13 760
Specific classifier 26188	2	5	5
TCM - Granulated herbal extracts (PhytoComm)	216	513	851
HCK - nutritional supplements (Hepart)	63	863	1881

From category *C* a total of 95 641 spectra from 30 708 batches for a total of 977 substances were taken into account for validation.

Validation results

The validation runs checked whether all substances/substance groups in the substance class *Narcotic substances, solid* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, the matching of all relevant spectra of various substances with the substances/substance groups in the substance class *Narcotic substances, solid* was checked and the correctness of the results was evaluated. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	6941	10	1 616 784
Type B	0	3537	112	3 640 841
Type C	0	995	8	1 433 612

All substances/substance groups in the substance class *Narcotic substances, solid* can be clearly distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.000 00 % (> 99.999 39 %)	99.785 64 % (> 99.735 38 %)
Type B	100.000 00 % (> 99.999 57 %)	97.441 93 % (> 97.331 26 %)
Type C	100.000 00 % (> 99.991 95 %)	95.980 27 % (> 91.346 41 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Amphetamine sulfate**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 20299-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Amphetamine sulfate; Amfetamini sulfas; Amphetamine sulphate

Special notes

When selecting the *Amphetamine sulfate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Amphetamine sulfate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Amphetamine sulf...	13G08-N11	21130	40	0814K-03361
Caelo	Amphetamine sulf...	13344603	21132	40	0814J-03361
Th. Geyer	Amphetamine sulf...	67002	21134	40	0814I-03361
Fagron	Amphetamine sulf...	15L01-B05-315291	22430	60	1316Q-03361
Th. Geyer	Amphetamine sulf...	067043	23459	40	2017C-03361

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Amphetamine sulf. . .	16153902	23461	40	2017A-03361
Fagron	Amphetamine sulf. . .	13G08-N11	21130SI	60	0814K-03361
Caelo	Amphetamine sulf. . .	13344603	21132SI	60	0814J-03361
Th. Geyer	Amphetamine sulf. . .	67002	21134SI	60	0814I-03361
Fagron	Amphetamine sulf. . .	15L01-B05-315291	22430SI	60	1316Q-03361

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 500 spectra of 10 reference samples from the substance/substance group *Amphetamine sulfate*. These samples are listed above in the [calibration samples](#) section.
- 107 749 spectra from a total of 1896 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 240 spectra of 6 reference samples from the substance/substance group *Amphetamine sulfate*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Amphetamine sulfate	16A29-B03-324577	23460	40
Caelo	Amphetamine sulfate	18186907	24646	40
Th. Geyer	Amphetamine sulfate	067043	23459SI	40
Fagron	Amphetamine sulfate	16A29-B03-324577	23460SI	40
Caelo	Amphetamine sulfate	16153902	23461SI	40
Caelo	Amphetamine sulfate	18186907	24646SI	40

- 242 726 spectra from a total of 3903 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 210 spectra from 84 *Apo-Ident* customers from 56 batches from the substance/substance group *Amphetamine sulfate*.

Supplier	Substance	Batch	Spectra
Th. Geyer	Amphetamine sulfate	67002	45
Sanacorp	Amphetamine sulfate	12/24-N03	1
Fagron	Amphetamine sulfate	11/01-N01	1
Fagron	Amphetamine sulfate	12E11-N06	2
Fagron	Amphetamine sulfate	12B06-N13	2
Caelo	Amphetamine sulfate	12I24-N03	1
Th. Geyer	Amphetamine sulfate	21082013B	1
Phönix	Amphetamine sulfate	13B11N04	1
Caelo	Amphetamine sulfate	12B06-N013	1
Fagron	Amphetamine sulfate	11D04-N13	1
Fagron	Amphetamine sulfate	13B11-N04	1
Fagron	Amphetamine sulfate	12/24N-03	1
Gehe	Amphetamine sulfate	13B11-N04	1
Fagron	Amphetamine sulfate	13G08-N11	3
Caelo	Amphetamine sulfate	13344609	1
Caelo	Amphetamine sulfate	1ot067002	1
Gatt Koller/Herba	Amphetamine sulfate	942/03130814	1
Fagron	Amphetamine sulfate	13J18-B40-290538	1
Fagron	Amphetamine sulfate	13D09-N03	1
Fagron	Amphetamine sulfate	13J18-B40	1
Fagron	Amphetamine sulfate	78012	1
Fagron	Amphetamine sulfate	13J18B40	1
Fagron	Amphetamine sulfate	88689	2
Th.Geyer über Phoenix	Amphetamine sulfate	67002	2
Fagron	Amphetamine sulfate	14E19-B09	2
Hedinger	Amphetamine sulfate	67002	1
Fagron	Amphetamine sulfate	14E19-B09-294572	1
Fagron	Amphetamine sulfate	14E19-B09-304502	2
Caelo	Amphetamine sulfate	15193204	1
Th. Geyer	Amphetamine sulfate	PH10051	1
Th. Geyer	Amphetamine sulfate	67043	38
Th. Geyer	Amphetamine sulfate	L0T067043	1
Anzag/Th.Gyer GmbH	Amphetamine sulfate	67043	1
Caelo	Amphetamine sulfate	15L01-B05	2
Fagron	Amphetamine sulfate	15L01-B05-315291	3
Fagron	Amphetamine sulfate	15L01-B05	1
Fagron	Amphetamine sulfate	16A29-B03	1
Fagron	Amphetamine sulfate	2213Q-03361	1
Caelo	Amphetamine sulfate	67043	15
Fagron	Amphetamine sulfate	16A29B03318744	1
Fagron	Amphetamine sulfate	16A29-B03-318743	6
Fagron	Amphetamine sulfate	67043	1
Caelo	Amphetamine sulfate	16A29-B03	7
Fagron	Amphetamine sulfate	15101-b05-315292	1
Caelo	Amphetamine sulfate	16153903	4
Th. Geyer	Amphetamine sulfate	6703	1
Fagron	Amphetamine sulfate	15L01-B05-315290	1
Fagron	Amphetamine sulfate	96200	1
Fagron	Amphetamine sulfate	324577	2
Euro OTC	Amphetamine sulfate	67043	5
Fagron	Amphetamine sulfate	315291	1
Noweda	Amphetamine sulfate	16A29-B03-318744	2
BLS/Noweda	Amphetamine sulfate	67043	1
Caelo	Amphetamine sulfate	16153909	1
Fagron	Amphetamine sulfate	16A29-B03-324578	1
Caelo	Amphetamine sulfate	16153902	1
Alkoholvertrieb Süd	Amphetamine sulfate	67043	2
Th. Geyer	Amphetamine sulfate	20170808-03361	1
Fagron	Amphetamine sulfate	16A29-B03-328317	3
Kehr	Amphetamine sulfate	67043	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Amphetamine sulfate	67043	1
Phönix	Amphetamine sulfate	16a29-b03-328317	1
Noweda	Amphetamine sulfate	16A29-B03-324578	1
Caelo	Amphetamine sulfate	1807171	1
Caelo	Amphetamine sulfate	16A29-B03-318744	1
Phönix	Amphetamine sulfate	16A29-B03-328317	1
Fagron	Amphetamine sulfate	16K08-B04-336361	1
Caelo	Amphetamine sulfate	17K03-B01	1
L30	Amphetamine sulfate	16A29-B03-328317	1
Fagron	Amphetamine sulfate	17K03-B01-344122	1
Th. Geyer	Amphetamine sulfate	70325	9
Th. Geyer	Amphetamine sulfate	57043	1
Caelo	Amphetamine sulfate	16143001	1

- 95 431 spectra from 646 *Apo-Ident* customers from a total of 30 139 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Amphetamine sulfate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Amphetamine sulfate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	500	0	107 749
Type B	0	235	5	242 726
Type C	0	208	2	95 431

The substance/substance group *Amphetamine sulfate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 98.8000 %)
Type B	100.0000 % (> 99.9935 %)	97.9167 % (> 96.6667 %)
Type C	100.0000 % (> 99.8791 %)	99.0476 % (> 97.6190 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
23459	23459	0.00	36.32
23461	23461	0.00	36.65
21130	21130	0.00	35.29
21132	21132	0.00	35.23
21134	21134	0.00	34.85
22430	22430	0.00	37.81

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Barbital sodium**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 20309-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Barbital sodium; Barbitalum natriicum

Special notes

When selecting the *Barbital sodium* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Barbital sodium*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Barbital sodium	12F27-N02	21136	40	0814I-07775
Fagron	Barbital sodium	14K14-B04-303353	22429	60	1516I-07775
Fagron	Barbital sodium	12F27-N02	21136SI	60	0814I-07775
Fagron	Barbital sodium	14K14-B04-303353	22429SI	60	1516I-07775

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 220 spectra of 4 reference samples from the substance/substance group *Barbital sodium*. These samples are listed above in the [calibration samples](#) section.
- 108 029 spectra from a total of 1900 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 80 spectra of 2 reference samples from the substance/substance group *Barbital sodium*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Barbital sodium	15L07-B02-324539	23463	40
Fagron	Barbital sodium	15L07-B02-324539	23463SI	40

- 242 886 spectra from a total of 3906 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 1 spectra from 1 *Apo-Ident* customers from 1 batches from the substance/substance group *Barbital sodium*.

Supplier	Substance	Batch	Spectra
Fagron	Barbital sodium	15L07-B02	1

- 95 640 spectra from 646 *Apo-Ident* customers from a total of 30 194 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Barbital sodium* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Barbital sodium* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	220	0	108 029
Type B	0	80	0	242 886
Type C	0	1	0	95 640

The substance/substance group *Barbital sodium* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 97.2727 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 92.5000 %)
Type C	100.0000 % (> 99.8798 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21136	21136	0.00	81.21
22429	22429	0.00	77.13

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Cocaine hydrochloride**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 21140-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Cocaine hydrochloride; Cocainum hydrochloricum

Special notes

When selecting the *Cocaine hydrochloride* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Cocaine hydrochloride*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cocaine hydrochl...	13259305	21140	40	0814J-01278
Fagron	Cocaine hydrochl...	13I17-N02	21142	40	0814I-01278
Th. Geyer	Cocaine hydrochl...	60145	21144	40	0814K-01278
Caelo	Cocaine hydrochl...	15224802	22406	60	1316Q-01278
Caelo	Cocaine hydrochl...	15224802	22411	60	1316Q-01278

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Cocaine hydrochl...	17003902	23470	40	1917S-01278
Fagron	Cocaine hydrochl...	16B11-B05-323386	23471	40	1917Q-1278
Caelo	Cocaine hydrochl...	13259305	21140SI	60	0814J-01278
Fagron	Cocaine hydrochl...	13I17-N02	21142SI	60	0814I-01278
Th. Geyer	Cocaine hydrochl...	60145	21144SI	60	0814K-01278
Caelo	Cocaine hydrochl...	15224802	22406SI	60	1316Q-01278
Caelo	Cocaine hydrochl...	15224802	22411SI	60	1316Q-01278

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 620 spectra of 12 reference samples from the substance/substance group *Cocaine hydrochloride*. These samples are listed above in the [calibration samples](#) section.
- 107 629 spectra from a total of 1896 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 400 spectra of 10 reference samples from the substance/substance group *Cocaine hydrochloride*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Cocaine hydrochloride	14H27-B03-304551	23472	40
Caelo	Cocaine hydrochloride	16142710	23473	40
Caelo	Cocaine hydrochloride	18266104	24648	40
Fagron	Cocaine hydrochloride	18J01-B07-353630	24649	40
Caelo	Cocaine hydrochloride	17003902	23470SI	40
Fagron	Cocaine hydrochloride	16B11-B05-323386	23471SI	40
Fagron	Cocaine hydrochloride	14H27-B03-304551	23472SI	40
Caelo	Cocaine hydrochloride	16142710	23473SI	40
Caelo	Cocaine hydrochloride	18266104	24648SI	40
Fagron	Cocaine hydrochloride	18J01-B07-353630	24649SI	40

- 242 566 spectra from a total of 3901 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 23 spectra from 11 *Apo-Ident* customers from 17 batches from the substance/substance group *Cocaine hydrochloride*.

Supplier	Substance	Batch	Spectra
Fagron	Cocaine hydrochloride	12G23-N13	1
Caelo	Cocaine hydrochloride	12329103	1
Fagron	Cocaine hydrochloride	13G03-N08	2
Caelo	Cocaine hydrochloride	13112102	1
Caelo	Cocaine hydrochloride	14069209	1
Caelo	Cocaine hydrochloride	13259306	1
Fagron	Cocaine hydrochloride	14H27-B03	1
Fagron	Cocaine hydrochloride	15E12-B03	1
Macfarlan Smith/Noweda	Cocaine hydrochloride	15224802	2
Synopharm	Cocaine hydrochloride	16B11-B05	1
Caelo	Cocaine hydrochloride	17003902	2
Fagron	Cocaine hydrochloride	16B11-B05-323385	1
Fagron	Cocaine hydrochloride	15-00638	1
Caelo	Cocaine hydrochloride	16051402	1
Gatt-Koller	Cocaine hydrochloride	20172282	1
Caesar & Loretz GmbH/Ph...	Cocaine hydrochloride	18266104	1
Caelo	Cocaine hydrochloride	18266104	2
Caelo	Cocaine hydrochloride	18J01-B07-353630	2

- 95 618 spectra from 646 *Apo-Ident* customers from a total of 30 178 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Cocaine hydrochloride* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Cocaine hydrochloride* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	620	0	107 629
Type B	0	400	0	242 566
Type C	0	21	2	95 618

The substance/substance group *Cocaine hydrochloride* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 99.0323 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 98.5000 %)
Type C	100.0000 % (> 99.8791 %)	91.3043 % (> 78.2609 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21140	21140	0.00	33.51
21142	21142	0.00	34.62
21144	21144	0.00	34.60
22406	22406	0.00	32.71
23470	23470	0.00	30.57
22411	22411	0.00	33.98
23471	23471	0.00	33.77

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group	Codeine phosphate hemihydrate
Substance class	Narcotic substances, solid
Report date	06/12/2019
Report number	20302-2019-12-06
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Codeine phosphate hemihydrate; Codeini phosphas hemihydricus

Special notes

When selecting the *Codeine phosphate hemihydrate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5	<i>European Pharmacopoeia 9th Edition, Basic Version 2017</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Appendix F	Addition to the second-stage models

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Codeine phosphate hemihydrate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Codeine phosphat...	13C11-N10	21148	40	0814I-01279
Th. Geyer	Codeine phosphat...	67038	21150	40	0814K-01279
Caelo	Codeine phosphat...	13071504	21152	40	0814J-01279
Caelo	Codeine phosphat...	15350703	22407	60	1316R-01279
Fagron	Codeine phosphat...	15F09-B02-310520	22425	60	1316Q-01279

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Th. Geyer	Codeine phosphat...	067054	23466	40	2017D-01279
Fagron	Codeine phosphat...	16E17-B05-324109	23469	40	2017A-01279
Fagron	Codeine phosphat...	13C11-N10	21148SI	60	0814I-01279
Th. Geyer	Codeine phosphat...	67038	21150SI	60	0814K-01279
Caelo	Codeine phosphat...	13071504	21152SI	60	0814J-01279
Caelo	Codeine phosphat...	15350703	22407SI	60	1316R-01279
Fagron	Codeine phosphat...	15F09-B02-310520	22425SI	60	1316Q-01279

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 620 spectra of 12 reference samples from the substance/substance group *Codeine phosphate hemihydrate*. These samples are listed above in the [calibration samples](#) section.
- 107 629 spectra from a total of 1895 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 320 spectra of 8 reference samples from the substance/substance group *Codeine phosphate hemihydrate*.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Codeine phosphate hemihydrate	16349702	23467	40
Caelo	Codeine phosphate hemihydrate	16169807	23468	40
Caelo	Codeine phosphate hemihydrate	18187906	24650	40
Th. Geyer	Codeine phosphate hemihydrate	067054	23466SI	40
Caelo	Codeine phosphate hemihydrate	16349702	23467SI	40
Caelo	Codeine phosphate hemihydrate	16169807	23468SI	40
Fagron	Codeine phosphate hemihydrate	16E17-B05-324109	23469SI	40
Caelo	Codeine phosphate hemihydrate	18187906	24650SI	40

- 242 646 spectra from a total of 3902 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 84 spectra from 21 *Apo-Ident* customers from 34 batches from the substance/substance group *Codeine phosphate hemihydrate*.

Supplier	Substance	Batch	Spectra
Fagron	Codeine phosphate hemihydrate	11/19-N10	1
Fagron	Codeine phosphate hemihydrate	12C05-N01	1
Merck	Codeine phosphate hemihydrate	120765	1
Fagron	Codeine phosphate hemihydrate	13C11-N10	1
Fagron	Codeine phosphate hemihydrate	13E21-N04	1
Th. Geyer	Codeine phosphate hemihydrate	67038	1
Fagron	Codeine phosphate hemihydrate	13/27-N06	2
Fagron	Codeine phosphate hemihydrate	913I-01279	1
Caelo	Codeine phosphate hemihydrate	13I27-N06	1
Caelo	Codeine phosphate hemihydrate	1427910	1
Caelo	Codeine phosphate hemihydrate	15219805	2
Fagron	Codeine phosphate hemihydrate	15B17-B99-304623	1
Fagron	Codeine phosphate hemihydrate	15F09-B02-310520	2
Th. Geyer	Codeine phosphate hemihydrate	67054	31
Caelo	Codeine phosphate hemihydrate	15350703	2
Fagron	Codeine phosphate hemihydrate	1615I-01279	1
Caelo	Codeine phosphate hemihydrate	14275910	1
Gehe	Codeine phosphate hemihydrate	15350710	1
Caelo	Codeine phosphate hemihydrate	15350710	2
Fagron	Codeine phosphate hemihydrate	15F09-B02	1
Fagron	Codeine phosphate hemihydrate	16E17-B05-324109	4
Euro/ Sanacorp	Codeine phosphate hemihydrate	15F09-B02-315564	2
Caelo	Codeine phosphate hemihydrate	16169805	3
Caelo	Codeine phosphate hemihydrate	16349706	2
Caelo	Codeine phosphate hemihydrate	16349704	2
Fagron	Codeine phosphate hemihydrate	15F09-B02-314033	1
Caelo	Codeine phosphate hemihydrate	16349708	2
Caelo	Codeine phosphate hemihydrate	17160503	1
Euro OTC	Codeine phosphate hemihydrate	CP04236	1
Caelo	Codeine phosphate hemihydrate	130715	1
Caelo	Codeine phosphate hemihydrate	18187903	1
Caesar & Loretz GmbH	Codeine phosphate hemihydrate	18187903	1
Caesar & Loretz GmbH	Codeine phosphate hemihydrate	18187906	2
Caelo	Codeine phosphate hemihydrate	18187902	1
Caesar & Loretz GmbH	Codeine phosphate hemihydrate	18187902	1
Caelo	Codeine phosphate hemihydrate	18187906	1
Caelo	Codeine phosphate hemihydrate	18187960	2
Caesar & Loretz GmbH	Codeine phosphate hemihydrate	18187901	1

- 95 557 spectra from 646 *Apo-Ident* customers from a total of 30 161 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Codeine phosphate hemihydrate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Codeine phosphate*

hemihydrate and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	620	0	107 629
Type B	0	319	1	242 646
Type C	0	84	0	95 557

The substance/substance group *Codeine phosphate hemihydrate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 99.0323 %)
Type B	100.0000 % (> 99.9935 %)	99.6875 % (> 98.7500 %)
Type C	100.0000 % (> 99.8791 %)	100.0000 % (> 92.8571 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21150	21150	0.00	51.09
21152	21152	0.00	54.13
22407	22407	0.00	52.70
21148	21148	0.00	52.42
22425	22425	0.00	52.17
23466	23466	0.00	50.08
23469	23469	0.00	55.02

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Dexamphetamine sulfate**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 21154-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Dexamphetamine sulfate; Dexamphetamine sulphate; Dexamphetaminum sulfuricum

Special notes

When selecting the *Dexamphetamine sulfate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dexamphetamine sulfate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Dexamphetamine s...	13I23-N01	21154	40	0814M-03582
Fagron	Dexamphetamine s...	15H06-B04-313862	22426	60	1316Q-03582
Fagron	Dexamphetamine s...	16D20-B01-321049	23474	40	2017B-03582
Fagron	Dexamphetamine s...	13I23-N01	21154SI	60	0814M-03582
Fagron	Dexamphetamine s...	15H06-B04-313862	22426SI	60	1316Q-03582

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 260 spectra of 5 reference samples from the substance/substance group *Dexamphetamine - sulfate*. These samples are listed above in the [calibration samples](#) section.
- 107 989 spectra from a total of 1899 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 120 spectra of 3 reference samples from the substance/substance group *Dexamphetamine - sulfate*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Dexamphetamine sulfate	16F09-B02-328352	23475	40
Fagron	Dexamphetamine sulfate	16D20-B01-321049	23474SI	40
Fagron	Dexamphetamine sulfate	16F09-B02-328352	23475SI	40

- 242 846 spectra from a total of 3905 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 8 spectra from 6 *Apo-Ident* customers from 6 batches from the substance/substance group *Dexamphetamine sulfate*.

Supplier	Substance	Batch	Spectra
Fargon/Gehe	Dexamphetamine sulfate	15H06-B04	1
Fargon	Dexamphetamine sulfate	11F29-N01	1
Caelo	Dexamphetamine sulfate	C9H13N14	1
Fargon	Dexamphetamine sulfate	15H06-B04-313861	1
Fargon	Dexamphetamine sulfate	17B09-B02-339689	3
Fargon	Dexamphetamine sulfate	17B09-B02-339690	1

- 95 633 spectra from 646 *Apo-Ident* customers from a total of 30 189 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dexamphetamine sulfate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dexamphetamine sulfate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	260	0	107 989
Type B	0	120	0	242 846
Type C	0	8	0	95 633

The substance/substance group *Dexamphetamine sulfate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 97.6923 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.8792 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21154	21154	0.00	21.04
22426	22426	0.00	24.49
23474	23474	0.00	23.33

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group	Dihydrocodeine hydrogen tartrate
Substance class	Narcotic substances, solid
Report date	06/12/2019
Report number	20303-2019-12-06
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Dihydrocodeine hydrogen tartrate; Dihydrocodeini hydrogenotartras

Special notes

When selecting the *Dihydrocodeine hydrogen tartrate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5	<i>European Pharmacopoeia 9th Edition, Basic Version 2017</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Appendix F	Addition to the second-stage models

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Dihydrocodeine hydrogen tartrate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Dihydrocodeine h...	13267601	21156	40	0914A-03065
Caelo	Dihydrocodeine h...	12142506	21158	40	0914B-03065
Fagron	Dihydrocodeine h...	13H19-N03	21160	40	0914C-03065
Th. Geyer	Dihydrocodeine h...	67035	21162	40	0914D-03065
Caelo	Dihydrocodeine h...	15295901	22418	60	1316Q-03065

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Dihydrocodeine h...	14K26-B04-306298	22422	60	1516I-03065
Th. Geyer	Dihydrocodeine h...	067052	23476	40	2017D-03065
Fagron	Dihydrocodeine h...	15F18-B07-323388	23479	40	2017A-03065
Caelo	Dihydrocodeine h...	13267601	21156SI	60	0914A-03065
Caelo	Dihydrocodeine h...	12142506	21158SI	60	0914B-03065
Fagron	Dihydrocodeine h...	13H19-N03	21160SI	60	0914C-03065
Th. Geyer	Dihydrocodeine h...	67035	21162SI	60	0914D-03065
Caelo	Dihydrocodeine h...	15295901	22418SI	60	1316Q-03065

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 660 spectra of 13 reference samples from the substance/substance group *Dihydrocodeine hydrogen tartrate*. These samples are listed above in the [calibration samples](#) section.
- 107 589 spectra from a total of 1894 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 380 spectra of 9 reference samples from the substance/substance group *Dihydrocodeine hydrogen tartrate*.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Dihydrocodeine hydrogen tart...	16106209	23477	40
Caelo	Dihydrocodeine hydrogen tart...	15295910	23478	40
Caelo	Dihydrocodeine hydrogen tart...	18115205	24651	40
Fagron	Dihydrocodeine hydrogen tart...	14K26-B04-306298	22422SI	60
Th. Geyer	Dihydrocodeine hydrogen tart...	067052	23476SI	40
Caelo	Dihydrocodeine hydrogen tart...	16106209	23477SI	40
Caelo	Dihydrocodeine hydrogen tart...	15295910	23478SI	40
Fagron	Dihydrocodeine hydrogen tart...	15F18-B07-323388	23479SI	40
Caelo	Dihydrocodeine hydrogen tart...	18115205	24651SI	40

- 242 586 spectra from a total of 3901 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 69 spectra from 12 *Apo-Ident* customers from 34 batches from the substance/substance group *Dihydrocodeine hydrogen tartrate*.

Supplier	Substance	Batch	Spectra
Th. Geyer	Dihydrocodeine hydrogen tart...	59819	1
Fagron	Dihydrocodeine hydrogen tart...	12G05-N01	1
Caelo	Dihydrocodeine hydrogen tart...	12142503	2
Fagron	Dihydrocodeine hydrogen tart...	13H19-N08	3
Th. Geyer	Dihydrocodeine hydrogen tart...	67035	2
Caelo	Dihydrocodeine hydrogen tart...	12142508	2
Caelo	Dihydrocodeine hydrogen tart...	13267603	1
Caelo	Dihydrocodeine hydrogen tart...	13267604	2
Caelo	Dihydrocodeine hydrogen tart...	13267608	1
Fagron	Dihydrocodeine hydrogen tart...	14D08-B06	1
Caelo	Dihydrocodeine hydrogen tart...	15295903	2
Fagron	Dihydrocodeine hydrogen tart...	14K26-B04	1
Caelo	Dihydrocodeine hydrogen tart...	15295907	1
Caelo	Dihydrocodeine hydrogen tart...	15295905	4
Fagron	Dihydrocodeine hydrogen tart...	14K26-B04-310974	1
Th. Geyer	Dihydrocodeine hydrogen tart...	67052	10
Fagron	Dihydrocodeine hydrogen tart...	14K26-B04-317257	1
Caelo	Dihydrocodeine hydrogen tart...	16106202	5
Fagron	Dihydrocodeine hydrogen tart...	15F18-B07	1
Caelo	Dihydrocodeine hydrogen tart...	15F18-B07-322054	1
Fagron	Dihydrocodeine hydrogen tart...	15F18-B07-323388	1
Hedinger	Dihydrocodeine hydrogen tart...	70162	4
Fagron	Dihydrocodeine hydrogen tart...	15F18B07	2
Caelo	Dihydrocodeine hydrogen tart...	16106203	4
Th. Geyer	Dihydrocodeine hydrogen tart...	20170808-03065	1
Fagron	Dihydrocodeine hydrogen tart...	16-00326	1
Fagron	Dihydrocodeine hydrogen tart...	16I27-B03-332122	1
Caelo	Dihydrocodeine hydrogen tart...	17097502	2
Fagron	Dihydrocodeine hydrogen tart...	16I27-B03	1
Fagron	Dihydrocodeine hydrogen tart...	17G04-B02-342861	1
Caelo	Dihydrocodeine hydrogen tart...	17097508	2
Fagron	Dihydrocodeine hydrogen tart...	17G04-B02	1
Caelo	Dihydrocodeine hydrogen tart...	18115202	3
Caelo	Dihydrocodeine hydrogen tart...	18115205	2

- 95 572 spectra from 646 *Apo-Ident* customers from a total of 30 161 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Dihydrocodeine hydrogen tartrate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Dihydrocodeine hydrogen tartrate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	660	0	107 589
Type B	0	378	2	242 586
Type C	0	69	0	95 572

The substance/substance group *Dihydrocodeine hydrogen tartrate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 99.0909 %)
Type B	100.0000 % (> 99.9935 %)	99.4737 % (> 98.6842 %)
Type C	100.0000 % (> 99.8791 %)	100.0000 % (> 91.3043 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21156	21156	0.00	31.41
21158	21158	0.00	31.02
21160	21160	0.00	30.76
21162	21162	0.00	31.34
22418	22418	0.00	31.65
22422	22422	0.00	30.98
23476	23476	0.00	30.79
23479	23479	0.00	29.98

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Fentanyl citrate**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 21164-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Fentanyl citrate; Fentanili citras

Special notes

When selecting the *Fentanyl citrate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Fentanyl citrate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Fentanyl citrate	13G30-N06	21164	40	0814M-07678
Fagron	Fentanyl citrate	15F18-B01-525098	22431	60	1316R-03686
Fagron	Fentanyl citrate	15F18-B01-525166	22432	60	1316Q-03686
Fagron	Fentanyl citrate	13G30-N06	21164SI	60	0814M-07678
Fagron	Fentanyl citrate	15F18-B01-525098	22431SI	60	1316R-03686
Fagron	Fentanyl citrate	15F18-B01-525166	22432SI	60	1316Q-03686

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 340 spectra of 6 reference samples from the substance/substance group *Fentanyl citrate*. These samples are listed above in the [calibration samples](#) section.
- 107 909 spectra from a total of 1899 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 80 spectra of 2 reference samples from the substance/substance group *Fentanyl citrate*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Fentanyl citrate	16E26-B01-525360	23480	40
Fagron	Fentanyl citrate	16E26-B01-525360	23480SI	40

- 242 886 spectra from a total of 3906 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 2 spectra from 2 *Apo-Ident* customers from 2 batches from the substance/substance group *Fentanyl citrate*.

Supplier	Substance	Batch	Spectra
Fagron	Fentanyl citrate	20100911-1	1
Fagron	Fentanyl citrate	10H20-N02	1

- 95 639 spectra from 646 *Apo-Ident* customers from a total of 30 193 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#),

the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Fentanyl citrate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Fentanyl citrate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	340	0	107 909
Type B	0	80	0	242 886
Type C	0	1	1	95 639

The substance/substance group *Fentanyl citrate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 98.2353 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 92.5000 %)
Type C	100.0000 % (> 99.8794 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
22431	22431	0.00	54.13
22432	22432	0.00	53.99
21164	21164	0.00	44.95

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Hydromorphone hydrochloride**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 20296-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Hydromorphone hydrochloride; Hydromorphinonum hydrochloricum

Special notes

When selecting the *Hydromorphone hydrochloride* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Hydromorphone hydrochloride*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Hydromorphone hy...	13267702	21166	40	0914A-01334
Fagron	Hydromorphone hy...	13E16-N02	21168	40	0914B-01334
Caelo	Hydromorphone hy...	15393602	22404	60	1516J-01334
Fagron	Hydromorphone hy...	13I25-N05	22427	60	1516I-01334
Fagron	Hydromorphone hy...	15L08-B08-316877	23482	40	1917Q-01334

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Hydromorphone hy...	13267702	21166SI	60	0914A-01334
Fagron	Hydromorphone hy...	13E16-N02	21168SI	60	0914B-01334
Caelo	Hydromorphone hy...	15393602	22404SI	60	1516J-01334
Fagron	Hydromorphone hy...	13I25-N05	22427SI	60	1516I-01334

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 480 spectra of 9 reference samples from the substance/substance group *Hydromorphone hydrochloride*. These samples are listed above in the [calibration samples](#) section.
- 107 769 spectra from a total of 1897 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 280 spectra of 7 reference samples from the substance/substance group *Hydromorphone hydrochloride*.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Hydromorphone hydrochloride	16123804	23481	40
Fagron	Hydromorphone hydrochloride	17F01-B01-337510	24652	40
Caelo	Hydromorphone hydrochloride	18316002	24653	40
Caelo	Hydromorphone hydrochloride	16123804	23481SI	40
Fagron	Hydromorphone hydrochloride	15L08-B08-316877	23482SI	40
Fagron	Hydromorphone hydrochloride	17F01-B01-337510	24652SI	40
Caelo	Hydromorphone hydrochloride	18316002	24653SI	40

- 242 686 spectra from a total of 3903 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 11 spectra from 7 *Apo-Ident* customers from 8 batches from the substance/substance group *Hydromorphone hydrochloride*.

Supplier	Substance	Batch	Spectra
Fagron	Hydromorphone hydrochloride	12G05-N02	1
Caelo	Hydromorphone hydrochloride	15393605	1
Fagron	Hydromorphone hydrochloride	13I27-N05	1
Caelo	Hydromorphone hydrochloride	16123805	2
Caelo	Hydromorphone hydrochloride	14210405	1
Caelo	Hydromorphone hydrochloride	16273604	1
Fagron	Hydromorphone hydrochloride	17F0-B01-337510	2
Caelo	Hydromorphone hydrochloride	18316002	2

- 95 630 spectra from 646 *Apo-Ident* customers from a total of 30 187 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Hydromorphone hydrochloride* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Hydromorphone hydrochloride* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	480	0	107 769
Type B	0	280	0	242 686
Type C	0	11	0	95 630

The substance/substance group *Hydromorphone hydrochloride* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 98.7500 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 97.8571 %)
Type C	100.0000 % (> 99.8792 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21166	21166	0.00	51.50
21168	21168	0.00	51.75
22404	22404	0.00	56.08
22427	22427	0.00	52.00
23482	23482	0.00	52.01

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group	Methadone hydrochloride
Substance class	Narcotic substances, solid
Report date	06/12/2019
Report number	20306-2019-12-06
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Methadone hydrochloride; Methadoni hydrochloridum

Special notes

When selecting the *Methadone hydrochloride* substance/substance group, the following information is displayed to the user:

- Advice for the identification of methadone hydrochloride: Levomethadone hydrochloride is a similar substance, which is not included in the database. An additional analysis of identity may be required, e.g. a melting-point measurement!

Applicable documents

978-3-7692-6679-5	<i>European Pharmacopoeia 9th Edition, Basic Version 2017</i> [3]
Komm2.2.40	<i>Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident</i> [4]
AA004	<i>Erstellung und Validierung eines IdentModul-Updates</i>
Appendix F	Addition to the second-stage models

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Methadone hydrochloride*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Methadone hydroc...	12344818	21170	40	0914E-02949
Fagron	Methadone hydroc...	13J24-B06	21172	40	0914F-02949
Th. Geyer	Methadone hydroc...	67033	21174	40	0914G-02949

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Th. Geyer	Methadone hydroc...	67047	21176	40	0914H-02949
Caelo	Methadone hydroc...	15244715	22417	60	1516A-02949
Fagron	Methadone hydroc...	15E29-B02-313098	22428	60	1516B-02949
Fagron	Methadone hydroc...	16F24-B02-323183	23485	40	1917Q-02949
Caelo	Methadone hydroc...	12344818	21170SI	60	0914E-02949
Fagron	Methadone hydroc...	13J24-B06	21172SI	60	0914F-02949
Th. Geyer	Methadone hydroc...	67033	21174SI	60	0914G-02949
Th. Geyer	Methadone hydroc...	67047	21176SI	60	0914H-02949
Caelo	Methadone hydroc...	15244715	22417SI	60	1516A-02949
Fagron	Methadone hydroc...	15E29-B02-313098	22428SI	60	1516B-02949

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 680 spectra of 13 reference samples from the substance/substance group *Methadone hydrochloride*. These samples are listed above in the [calibration samples](#) section.
- 107 569 spectra from a total of 1895 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 360 spectra of 9 reference samples from the substance/substance group *Methadone hydrochloride*.

Supplier	Substance	Batch	Sample ID	Spectra
Th. Geyer	Methadone hydrochloride	067055	23483	40
Caelo	Methadone hydrochloride	16135309	23484	40
Fagron	Methadone hydrochloride	18G02-B02-355990	24654	40
Caelo	Methadone hydrochloride	18132413	24655	40
Th. Geyer	Methadone hydrochloride	067055	23483SI	40
Caelo	Methadone hydrochloride	16135309	23484SI	40
Fagron	Methadone hydrochloride	16F24-B02-323183	23485SI	40
Fagron	Methadone hydrochloride	18G02-B02-355990	24654SI	40
Caelo	Methadone hydrochloride	18132413	24655SI	40

- 242 606 spectra from a total of 3902 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B*

section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 463 spectra from 150 *Apo-Ident* customers from 207 batches from the substance/substance group *Methadone hydrochloride*.

Supplier	Substance	Batch	Spectra
Fagron	Methadone hydrochloride	12H07-N07	1
Fagron	Methadone hydrochloride	20130912-6	1
Caelo	Methadone hydrochloride	10076407	1
Fagron	Methadone hydrochloride	12H07-N05	3
Caelo	Methadone hydrochloride	11227507	1
Caelo	Methadone hydrochloride	12171907	1
Caelo	Methadone hydrochloride	12171902	1
Fagron	Methadone hydrochloride	1213E-02949	1
Siegfried Ltd	Methadone hydrochloride	12B21-N03	1
Fagron	Methadone hydrochloride	9041306	1
Fagron	Methadone hydrochloride	171920-0009	1
Caelo	Methadone hydrochloride	12171914	1
Fagron	Methadone hydrochloride	12A05-No2	1
Fagron	Methadone hydrochloride	1225H008	1
Sanacorp	Methadone hydrochloride	12171913	1
Fagron	Methadone hydrochloride	13B06-N02	5
Fagron	Methadone hydrochloride	12B21N03	1
Caelo	Methadone hydrochloride	12344802	1
Sanacorp	Methadone hydrochloride	12344806	2
Fagron	Methadone hydrochloride	12b21-n03	1
Caelo	Methadone hydrochloride	12171906	2
Fagron	Methadone hydrochloride	13F17-N02	4
Fagron	Methadone hydrochloride	25091307	1
Caelo	Methadone hydrochloride	12H07-N05	1
Caelo	Methadone hydrochloride	12171918	1
Caelo	Methadone hydrochloride	13F17-N02	2
Caelo	Methadone hydrochloride	12344816	3
Caelo	Methadone hydrochloride	12344812	1
AHD	Methadone hydrochloride	12344812	1
Caelo	Methadone hydrochloride	17011406	1
Fagron	Methadone hydrochloride	11G11-N02	1
Caelo	Methadone hydrochloride	12344818	1
Caelo	Methadone hydrochloride	12344807	1
Caelo	Methadone hydrochloride	12344814	2
Fagron	Methadone hydrochloride	414022-1-P	1
Caelo	Methadone hydrochloride	12344819	1
Fagron	Methadone hydrochloride	1316H037	4
Th. Geyer	Methadone hydrochloride	67047	2
Fagron	Methadone hydrochloride	11D07-N09	1
Caelo	Methadone hydrochloride	13266502	2
Fagron	Methadone hydrochloride	13I24-B06-286908	1
Caelo	Methadone hydrochloride	13266504	2
Fagron	Methadone hydrochloride	13J24-B06-286907	1
Fagron	Methadone hydrochloride	25091406	1
Caelo	Methadone hydrochloride	1327H051	1
Caelo	Methadone hydrochloride	13266501	1
Fagron	Methadone hydrochloride	1327h051	1
VDL;87,00EUR;08.09.14	Methadone hydrochloride	12344821	1
Fagron	Methadone hydrochloride	13J24-B06-286908	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Methadone hydrochloride	13266505	4
Fagron	Methadone hydrochloride	14C24-B01-292857	4
Euro OTC	Methadone hydrochloride	14C24-B01	1
Bombastus	Methadone hydrochloride	14C24-B01-292858	6
Fagron	Methadone hydrochloride	4031517	1
Caelo	Methadone hydrochloride	13266511	2
Caelo	Methadone hydrochloride	13266515	1
Fagron	Methadone hydrochloride	14H15-B01-299381	1
Fagron	Methadone hydrochloride	29051502	1
Caelo	Methadone hydrochloride	13266518	3
Caelo	Methadone hydrochloride	13266514	1
Fagron	Methadone hydrochloride	14C24-B01	1
Fagron	Methadone hydrochloride	15B24-B02-304724	6
Fagron	Methadone hydrochloride	1409H043	1
Fagron	Methadone hydrochloride	8081501	1
Caelo	Methadone hydrochloride	13266521	1
Fagron	Methadone hydrochloride	1437H012	2
Hedinger/ Sanacorp	Methadone hydrochloride	15215505	3
Caelo	Methadone hydrochloride	15215501	1
Caelo	Methadone hydrochloride	15244704	3
Fagron	Methadone hydrochloride	14H15-B01	1
Fagron	Methadone hydrochloride	14H15-B01-299382	6
Fagron	Methadone hydrochloride	15244704	1
Caelo	Methadone hydrochloride	15244702	1
Caelo	Methadone hydrochloride	15244706	3
Fagron	Methadone hydrochloride	15E29-B02-309070	2
Caelo	Methadone hydrochloride	14C24-B01-292860	1
Caelo	Methadone hydrochloride	15244705	1
Fagron	Methadone hydrochloride	15E29B02309070	1
Caelo	Methadone hydrochloride	15244710	4
Faagron	Methadone hydrochloride	1437H011	1
Siegfried Ltd, Division...	Methadone hydrochloride	15244712	1
Caelo	Methadone hydrochloride	15244709	2
Fagron	Methadone hydrochloride	15I24-B07	1
Fagron	Methadone hydrochloride	15E29-B02-313031	1
Caelo	Methadone hydrochloride	15244701	1
Jenne	Methadone hydrochloride	15B24-B02	1
Fagron	Methadone hydrochloride	15E29-B02-313098	3
Caelo	Methadone hydrochloride	15244717	6
Caelo	Methadone hydrochloride	15244715	1
Fagron	Methadone hydrochloride	15B24-B02-313097	1
Caelo	Methadone hydrochloride	15244713	4
Caelo	Methadone hydrochloride	15215507	1
Fagron	Methadone hydrochloride	16B16-B02-318420	2
Fagron	Methadone hydrochloride	13J24B06286910	1
Fagron	Methadone hydrochloride	1514H042	3
Phoenix	Methadone hydrochloride	15244715	1
L30	Methadone hydrochloride	67055	6
Caelo	Methadone hydrochloride	16135302	1
Caelo	Methadone hydrochloride	15244718	1
Euro OTC	Methadone hydrochloride	16135304	5
Caelo	Methadone hydrochloride	15I24-B07-317266	1
Caelo	Methadone hydrochloride	15244720	1
Caelo	Methadone hydrochloride	16135308	3
Caelo	Methadone hydrochloride	16135307	2
Fagron	Methadone hydrochloride	1608H002	1
Caelo	Methadone hydrochloride	16135303	1
Caelo	Methadone hydrochloride	16135305	6
Fagron	Methadone hydrochloride	16F24-B02-323182	1
Caelo	Methadone hydrochloride	16135301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Methadone hydrochloride	1608H002	1
Fagron	Methadone hydrochloride	16F24-B02-323185	1
Fagron	Methadone hydrochloride	16K01-B04-329933	2
Fagron	Methadone hydrochloride	16F24-B02-326207	2
Fagron	Methadone hydrochloride	170427MH/1608H002	1
Dr. Wolff	Methadone hydrochloride	16F24-B02	1
Caelo	Methadone hydrochloride	28041725	1
Caelo	Methadone hydrochloride	16135310	4
Phönix	Methadone hydrochloride	16135309	1
Fagron	Methadone hydrochloride	16K01B04	1
Fiebig/	Methadone hydrochloride	67055	1
Fagron	Methadone hydrochloride	16F24-B02-323186	2
Th. Geyer	Methadone hydrochloride	67055	6
Caelo	Methadone hydrochloride	16F24-B02-323184	1
Caelo	Methadone hydrochloride	16135312	3
Fagron	Methadone hydrochloride	16B16-Bo2	1
Caelo	Methadone hydrochloride	16135311	4
BLS-Analytik/Sanacorp	Methadone hydrochloride	67055	3
Caelo	Methadone hydrochloride	28061702	1
Fagron	Methadone hydrochloride	16K01-B04-329934	1
Caelo	Methadone hydrochloride	16135318	2
Caelo	Methadone hydrochloride	16135315	1
Ebert-Jacobi	Methadone hydrochloride	16135306	1
Caelo	Methadone hydrochloride	16135317	1
Caelo	Methadone hydrochloride	16135309	1
Fagron	Methadone hydrochloride	16K01-B04	2
Euro OTC	Methadone hydrochloride	16135318	1
Caelo	Methadone hydrochloride	16135324	3
Gehe/Noweda	Methadone hydrochloride	16135319	1
Caelo	Methadone hydrochloride	16135326	5
Phönix	Methadone hydrochloride	67055	1
Caelo	Methadone hydrochloride	20170728/Ma2	1
Fagron	Methadone hydrochloride	16F24B02323184	2
Caelo	Methadone hydrochloride	16135325	1
Caelo	Methadone hydrochloride	16135327	9
Fagron	Methadone hydrochloride	16f24-b02-323184	1
Caelo	Methadone hydrochloride	16135319	1
Th. Geyer	Methadone hydrochloride	70245	30
Fargon	Methadone hydrochloride	16K01-B04-329937	1
Caelo	Methadone hydrochloride		1
Bombastus	Methadone hydrochloride	16135326	1
Fagron	Methadone hydrochloride	16K01-B04-329937	1
Caelo	Methadone hydrochloride	16135329	6
Caelo	Methadone hydrochloride	1613529	1
Fagron	Methadone hydrochloride	17C24-B05-334309	5
Caelo	Methadone hydrochloride	7091708	1
Caelo	Methadone hydrochloride	16135330	4
gehe	Methadone hydrochloride	16135316	2
Fiebig/BLS	Methadone hydrochloride	70245	1
Fagron	Methadone hydrochloride	16-k01-B04-329937	1
Caelo	Methadone hydrochloride	16135332	7
Caelo	Methadone hydrochloride	17C24-B05	1
Caelo	Methadone hydrochloride	16135316	1
Caelo	Methadone hydrochloride	17239504	3
Caelo	Methadone hydrochloride	16208016	1
Fagron	Methadone hydrochloride	17C24-B05	1
Fagron	Methadone hydrochloride	16F24-B02-329935	1
Caelo	Methadone hydrochloride	25111701	1
L30	Methadone hydrochloride	70245	1
Caelo	Methadone hydrochloride	172329504	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
BLS Analytik GmbH	Methadone hydrochloride	70245	7
Th. Geyer	Methadone hydrochloride	70246	1
Gehe/P.W.Beyvers	Methadone hydrochloride	70245	4
Caelo	Methadone hydrochloride	17239502	2
Fagron	Methadone hydrochloride	17I15-B02-339529	2
Caelo	Methadone hydrochloride	2912161	1
Caelo	Methadone hydrochloride	1729504	1
gehe	Methadone hydrochloride	17239502	1
Fagron	Methadone hydrochloride	17C24-B05-335711	7
Fagron	Methadone hydrochloride	17c24b05334309	1
Caelo	Methadone hydrochloride	17239506	5
Fagron	Methadone hydrochloride	17I15-B02-339521	2
BLS	Methadone hydrochloride	70245	3
Fagron	Methadone hydrochloride	1718H028	3
Fagro//Noweda	Methadone hydrochloride	15/24-B07-317267	1
Phönix	Methadone hydrochloride	70245	1
Caelo	Methadone hydrochloride	16135323	1
Caelo	Methadone hydrochloride	16135328	1
Caelo	Methadone hydrochloride	17239507	3
Caelo	Methadone hydrochloride	180219MC/17239506	1
Caelo	Methadone hydrochloride	23021802	1
Caelo	Methadone hydrochloride	17239503	3
Fagron	Methadone hydrochloride	17/15-B02-339521	2
Th. Geyer	Methadone hydrochloride		1
Fagron	Methadone hydrochloride	17C24-B05-334311	2
Caelo	Methadone hydrochloride	17I15-B02-339523	2
Caelo	Methadone hydrochloride	17239501	1
Caelo	Methadone hydrochloride	17239510	2
Fagro/Jenne	Methadone hydrochloride	17I15-B02-339525	2
Fagro/ Noweda	Methadone hydrochloride	17/15-B02-339527	1
Caelo	Methadone hydrochloride	17239512	7
Fagron	Methadone hydrochloride	17I15W02	1
Caelo	Methadone hydrochloride	17051801	1
Caelo	Methadone hydrochloride	1723	1
Fagron	Methadone hydrochloride	17I15-B02-339523	2
Caelo	Methadone hydrochloride	17239518	2
Caelo	Methadone hydrochloride	17239514	1
Fagron	Methadone hydrochloride	17I15-B02	2
Caelo	Methadone hydrochloride	1081801	1
Fagron	Methadone hydrochloride	17I15-B02-343153	1
Caelo	Methadone hydrochloride	171920-0004	1
Fagron	Methadone hydrochloride	17I15-B02-339521	1
Caelo	Methadone hydrochloride	18186601	3
Fagron	Methadone hydrochloride	17/15-B02-343153	1
Caelo	Methadone hydrochloride	1723509	1
Caelo	Methadone hydrochloride	18186604	1
Caelo	Methadone hydrochloride	3111806	1
Caelo	Methadone hydrochloride	18132403	2
Fagron	Methadone hydrochloride	18C01-B05-350571	3
Th. Geyer	Methadone hydrochloride	56(070245)	1
Fagron	Methadone hydrochloride	18G02-B02	1
Caelo	Methadone hydrochloride	19011913	1
Fagron GmbH & Co. KG	Methadone hydrochloride	18C01-B05-350571	1
Gatt-Koller	Methadone hydrochloride	1801702	1
Caelo	Methadone hydrochloride	31011901	1
Caelo	Methadone hydrochloride	18132405	2
Fagron	Methadone hydrochloride	18C01B05350571	1
Caelo	Methadone hydrochloride	18132408	3
Fagron	Methadone hydrochloride	18G02-B02-351285	1
Caelo	Methadone hydrochloride	18132412	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Methadone hydrochloride	350571	1
Caelo	Methadone hydrochloride	57(18132413)	1
Fagron	Methadone hydrochloride	18G02-B02-355990	2
Caelo	Methadone hydrochloride	18132413	3
Caelo	Methadone hydrochloride	58(18132417)	1
Caesar & Loretz GmbH	Methadone hydrochloride	18132413	1

- 95 178 spectra from 646 *Apo-Ident* customers from a total of 29 990 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Methadone hydrochloride* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Methadone hydrochloride* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	680	0	107 569
Type B	0	360	0	242 606
Type C	0	460	3	95 178

The substance/substance group *Methadone hydrochloride* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 99.1176 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 98.3333 %)
Type C	100.0000 % (> 99.8791 %)	99.3521 % (> 98.7041 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21170	21170	0.00	43.05
21172	21172	0.00	43.23
21174	21174	0.00	32.96
21176	21176	0.00	32.67
22417	22417	0.00	44.36
22428	22428	0.00	43.57
23485	23485	0.00	43.89

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Midazolam**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 21178-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Midazolam; Midazolamum

Special notes

When selecting the *Midazolam* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Midazolam*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Midazolam	13K05-N01	21178	40	0914E-07393
Fagron	Midazolam	14C12-B31-313222	22424	60	1516A-07393
Fagron	Midazolam	13K05-N01	21178SI	60	0914E-07393
Fagron	Midazolam	14C12-B31-313222	22424SI	60	1516A-07393

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 220 spectra of 4 reference samples from the substance/substance group *Midazolam*. These samples are listed above in the [calibration samples](#) section.
- 108 029 spectra from a total of 1900 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 150 spectra of 4 reference samples from the substance/substance group *Midazolam*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Midazolam	16J03-B09-326963	23486	40
Fagron	Midazolam	17J16-B01-342518	24656	40
Fagron	Midazolam	16J03-B09-326963	23486SI	30
Fagron	Midazolam	17J16-B01-342518	24656SI	40

- 242 816 spectra from a total of 3905 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 2 spectra from 1 *Apo-Ident* customers from 2 batches from the substance/substance group *Midazolam*.

Supplier	Substance	Batch	Spectra
Fagron	Midazolam	14C12-B31-313222	1
Fagron	Midazolam	16J03-B09-326963	1

- 95 639 spectra from 646 *Apo-Ident* customers from a total of 30 193 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are docu-

mented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Midazolam* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Midazolam* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	220	0	108 029
Type B	0	135	15	242 816
Type C	0	2	0	95 639

The substance/substance group *Midazolam* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 97.2727 %)
Type B	100.0000 % (> 99.9935 %)	90.0000 % (> 88.0000 %)
Type C	100.0000 % (> 99.8794 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21178	21178	0.00	58.31
22424	22424	0.00	59.07

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Midazolam hydrochloride**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 20307-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Midazolam hydrochloride; Midazolami hydrochloridum

Special notes

When selecting the *Midazolam hydrochloride* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Midazolam hydrochloride*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Midazolam hydroc...	12215811	21180	40	0914E-03538
Fagron	Midazolam hydroc...	13E22-N09	21182	40	0914F-03538
Caelo	Midazolam hydroc...	15224602	22405	60	1516B-03538
Fagron	Midazolam hydroc...	14L23-B04-300245	22423	60	1516A-03538
Fagron	Midazolam hydroc...	15H04-B03-314396	23489	40	1917R-03538

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Midazolam hydroc...	16C14-B12-332469	23925	40	0718E-03538
Caelo	Midazolam hydroc...	12215811	21180SI	60	0914E-03538
Fagron	Midazolam hydroc...	13E22-N09	21182SI	60	0914F-03538
Caelo	Midazolam hydroc...	15224602	22405SI	60	1516B-03538
Fagron	Midazolam hydroc...	14L23-B04-300245	22423SI	60	1516A-03538
Fagron	Midazolam hydroc...	16C14-B12-332469	23925SI	40	0718E-03538

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 560 spectra of 11 reference samples from the substance/substance group *Midazolam hydrochloride*. These samples are listed above in the [calibration samples](#) section.
- 107 689 spectra from a total of 1896 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 250 spectra of 7 reference samples from the substance/substance group *Midazolam hydrochloride*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Midazolam hydrochloride	15J07-B01-320717	23487	40
Caelo	Midazolam hydrochloride	15224611	23488	40
Fagron	Midazolam hydrochloride	18B28-B05-354856	24657	40
Fagron	Midazolam hydrochloride	15J07-B01-320717	23487SI	30
Caelo	Midazolam hydrochloride	15224611	23488SI	30
Fagron	Midazolam hydrochloride	15H04-B03-314396	23489SI	30
Fagron	Midazolam hydrochloride	18B28-B05-354856	24657SI	40

- 242 716 spectra from a total of 3903 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 36 spectra from 20 *Apo-Ident* customers from 21 batches from the substance/substance group *Midazolam hydrochloride*.

Supplier	Substance	Batch	Spectra
Caelo	Midazolam hydrochloride	11184705	1
Eurho/Phoenix	Midazolam hydrochloride	12L17-N10	1
Fagron	Midazolam hydrochloride	9D27-N01	2
Fagron	Midazolam hydrochloride	14H05-B08-297982	1
Caelo	Midazolam hydrochloride	14339009	1
Caelo	Midazolam hydrochloride	14339015	2
Fagron	Midazolam hydrochloride	11A31-N01	1
Caelo	Midazolam hydrochloride	15224603	1
Caelo	Midazolam hydrochloride	15224605	1
Caelo	Midazolam hydrochloride	15224607	1
Caelo	Midazolam hydrochloride	15224614	1
Euro OTC	Midazolam hydrochloride	15224607	1
Caelo	Midazolam hydrochloride	14339002	1
Fagron	Midazolam hydrochloride	15H04-B03-314396	6
Caelo	Midazolam hydrochloride	15224616	2
Fagron	Midazolam hydrochloride	400101	1
Caelo	Midazolam hydrochloride	16290605	3
Fagron	Midazolam hydrochloride	16C14-B12-332469	1
Caelo	Midazolam hydrochloride	17168003	3
Caelo	Midazolam hydrochloride	83600116001	1
Caelo/ Noweda	Midazolam hydrochloride	17168006	1
Sanacorp	Midazolam hydrochloride	16290605	1
Sanacorp	Midazolam hydrochloride	17168006	1
Fagron	Midazolam hydrochloride	18B28-B05-354856	1

- 95 605 spectra from 646 *Apo-Ident* customers from a total of 30 174 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Midazolam hydrochloride* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Midazolam hydrochloride* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	560	0	107 689
Type B	0	250	0	242 716
Type C	0	36	0	95 605

The substance/substance group *Midazolam hydrochloride* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 98.9286 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 97.6000 %)
Type C	100.0000 % (> 99.8791 %)	100.0000 % (> 83.3333 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
23489	23489	0.00	72.03
23925	23925	0.00	52.06
21180	21180	0.00	86.37
21182	21182	0.00	85.80
22405	22405	0.00	75.89
22423	22423	0.00	77.18

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Morphine hydrochloride**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 21190-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Morphine hydrochloride; Morphin hydrochloridum

Special notes

When selecting the *Morphine hydrochloride* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Morphine hydrochloride*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Morphine hydroch...	13125207	21190	40	0914E-01879
Caelo	Morphine hydroch...	13349505	21192	40	0914F-01879
Fagron	Morphine hydroch...	13G08-N04	21194	40	0914G-01879
Th. Geyer	Morphine hydroch...	67036	21196	40	0914H-01879
Caelo	Morphine hydroch...	15439404	22409	60	1416R-01879

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Fagron	Morphine hydroch...	15I09-B08-312694	22421	60	1416Q-01879
Caelo	Morphine hydroch...	16184605	23490	40	1917S-01879
Caelo	Morphine hydroch...	13125207	21190SI	60	0914E-01879
Caelo	Morphine hydroch...	13349505	21192SI	60	0914F-01879
Fagron	Morphine hydroch...	13G08-N04	21194SI	60	0914G-01879
Th. Geyer	Morphine hydroch...	67036	21196SI	60	0914H-01879
Caelo	Morphine hydroch...	15439404	22409SI	60	1416R-01879
Fagron	Morphine hydroch...	15I09-B08-312694	22421SI	60	1416Q-01879

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 680 spectra of 13 reference samples from the substance/substance group *Morphine hydrochloride*. These samples are listed above in the [calibration samples](#) section.
- 107 569 spectra from a total of 1895 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 360 spectra of 9 reference samples from the substance/substance group *Morphine hydrochloride*.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Morphine hydrochloride	15439411	23491	40
Fagron	Morphine hydrochloride	16D21-B02-322088	23492	40
Fagron	Morphine hydrochloride	17F22-B01-338170	24658	40
Caelo	Morphine hydrochloride	18188607	24659	40
Caelo	Morphine hydrochloride	16184605	23490SI	40
Caelo	Morphine hydrochloride	15439411	23491SI	40
Fagron	Morphine hydrochloride	16D21-B02-322088	23492SI	40
Fagron	Morphine hydrochloride	17F22-B01-338170	24658SI	40
Caelo	Morphine hydrochloride	18188607	24659SI	40

- 242 606 spectra from a total of 3902 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 83 spectra from 34 *Apo-Ident* customers from 55 batches from the substance/substance group *Morphine hydrochloride*.

Supplier	Substance	Batch	Spectra
Caelo	Morphine hydrochloride	12318101	1
Fagron	Morphine hydrochloride	11A27-N04	1
Th. Geyer	Morphine hydrochloride	67036	5
Fagron	Morphine hydrochloride	12E10-N04	2
Caelo	Morphine hydrochloride	12104403	1
Caelo	Morphine hydrochloride	13081310	1
Caelo	Morphine hydrochloride	12104409	1
Caelo	Morphine hydrochloride	13125206	1
Fagron	Morphine hydrochloride	12E10N04	1
Fagron	Morphine hydrochloride	13H22-N12	1
Fagron	Morphine hydrochloride	13D09-N04	1
Th. Geyer	Morphine hydrochloride	PH120087	1
Caelo	Morphine hydrochloride	14055084	1
Caelo	Morphine hydrochloride	14238604	1
Caelo	Morphine hydrochloride	14238617	1
Fagron	Morphine hydrochloride	15A12-B04-303137	1
Th. Geyer	Morphine hydrochloride	1415E-01879	1
Caelo	Morphine hydrochloride	15199705	1
Caelo	Morphine hydrochloride		2
Caelo	Morphine hydrochloride	15199711	3
Caesar & Loretz GmbH	Morphine hydrochloride	15199711 (2)	1
Caesar & Loretz GmbH	Morphine hydrochloride	15199711 (1)	1
Caesar & Loretz GmbH	Morphine hydrochloride	15199712	1
Caesar & Loretz GmbH	Morphine hydrochloride	15199711	1
Synop	Morphine hydrochloride	15I09-B08	1
Fagron	Morphine hydrochloride	15I09-B08-312693	1
Macfarlan Smith Limited. .	Morphine hydrochloride	15199711	1
Caelo	Morphine hydrochloride	15199709	1
Caelo	Morphine hydrochloride	1543904	1
Fagron	Morphine hydrochloride	15J27B01316889	1
Caelo	Morphine hydrochloride	14238619	1
Th. Geyer	Morphine hydrochloride	2615M-01879	1
Caelo	Morphine hydrochloride	15439409	1
Caelo	Morphine hydrochloride	15439406	3
Caelo	Morphine hydrochloride	13125208	1
Caelo	Morphine hydrochloride	15439405	2
Caelo	Morphine hydrochloride	15439402	1
Caelo	Morphine hydrochloride	15199713	1
Caelo	Morphine hydrochloride	16184602	2
Caelo	Morphine hydrochloride	14081701	1
Caelo	Morphine hydrochloride	16184605	1
Caelo	Morphine hydrochloride	16184610	2
Sanacorp	Morphine hydrochloride	16184612	2
Caelo	Morphine hydrochloride	16184609	2
Caelo	Morphine hydrochloride	16184613	6
Fagron	Morphine hydrochloride	16D21-B02-322087	1
Caelo	Morphine hydrochloride	15199707	1
Caelo	Morphine hydrochloride	16184614	3
Sanacorp	Morphine hydrochloride	16184605	1
Caelo	Morphine hydrochloride	16184607	1
Fagron	Morphine hydrochloride	17C24-B03-334287	1
Caelo	Morphine hydrochloride	18188603	3
Caelo	Morphine hydrochloride	15439411	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Morphine hydrochloride	18188601	1
Fagron	Morphine hydrochloride	17F22-B01-338170	2
Caelo	Morphine hydrochloride	118188607	2
Caelo	Morphine hydrochloride	18188614	1

- 95 558 spectra from 646 *Apo-Ident* customers from a total of 30 140 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Morphine hydrochloride* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Morphine hydrochloride* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	680	0	107 569
Type B	0	354	6	242 606
Type C	0	83	0	95 558

The substance/substance group *Morphine hydrochloride* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 99.1176 %)
Type B	100.0000 % (> 99.9935 %)	98.3333 % (> 97.5000 %)
Type C	100.0000 % (> 99.8791 %)	100.0000 % (> 92.7711 %)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21190	21190	0.00	21.30

continued on the next page

continued from previous page

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21192	21192	0.00	24.73
21194	21194	0.00	20.53
21196	21196	0.00	20.74
22409	22409	0.00	19.79
22421	22421	0.00	19.77
23490	23490	0.00	20.19

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Morphine sulfate**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 21184-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Morphine sulfate; Morphine sulphate; Morphini sulfas

Special notes

When selecting the *Morphine sulfate* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Morphine sulfate*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Morphine sulfate	12153212	21184	40	0914F-07776
Caelo	Morphine sulfate	11333103	21186	40	0914E-07776
Fagron	Morphine sulfate	13J07-N04	21188	40	0914G-07776
Caelo	Morphine sulfate	15255703	22408	60	1416Q-07776
Fagron	Morphine sulfate	14F02-B06-315405	23494	40	1917Q-07776
Caelo	Morphine sulfate	11333103	21186SI	31	0914E-07776
Fagron	Morphine sulfate	13J07-N04	21188SI	60	0914G-07776

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 311 spectra of 7 reference samples from the substance/substance group *Morphine sulfate*. These samples are listed above in the [calibration samples](#) section.
- 107 938 spectra from a total of 1897 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 349 spectra of 8 reference samples from the substance/substance group *Morphine sulfate*.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Morphine sulfate	15295802	23493	40
Caelo	Morphine sulfate	17331510	24660	40
Caelo	Morphine sulfate	12153212	21184SI	60
Caelo	Morphine sulfate	11333103	21186SI [†]	29
Caelo	Morphine sulfate	15255703	22408SI	60
Caelo	Morphine sulfate	15295802	23493SI	40
Fagron	Morphine sulfate	14F02-B06-315405	23494SI	40
Caelo	Morphine sulfate	17331510	24660SI	40

- 242 617 spectra from a total of 3901 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 2 spectra from 1 *Apo-Ident* customers from 1 batches from the substance/substance group *Morphine sulfate*.

[†]These spectra were recorded from material of the same package as the corresponding calibration spectra (*Type A*). Hence, these spectra are not independent and are not of *Type B*.

Supplier	Substance	Batch	Spectra
Caelo	Morphine sulfate	17331510	2

- 95 639 spectra from 646 *Apo-Ident* customers from a total of 30 194 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Morphine sulfate* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Morphine sulfate* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	301	10	107 938
Type B	0	266	83	242 617
Type C	0	2	0	95 639

The substance/substance group *Morphine sulfate* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate* (*specificity*) and the weighted *true positive rate* (*recognition rate*) are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	96.7846 % (> 95.8199 %)
Type B	100.0000 % (> 99.9935 %)	76.2178 % (> 75.3582 %)
Type C	100.0000 % (> 99.8794 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21184	21184	0.00	35.16
21186	21186	0.00	35.98
21188	21188	0.00	28.51
22408	22408	0.00	38.40
23494	23494	0.00	24.17

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group	Oxycodone hydrochloride
Substance class	Narcotic substances, solid
Report date	06/12/2019
Report number	20310-2019-12-06
Executing company	HiperScan GmbH Weißeritzstraße 3 01067 Dresden Germany

Relevant substance names

Oxycodone hydrochloride; Oxycodoni hydrochloridum

Special notes

When selecting the *Oxycodone hydrochloride* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5	European Pharmacopoeia 9th Edition, Basic Version 2017 [3]
Komm2.2.40	Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident [4]
AA004	Erstellung und Validierung eines IdentModul-Updates
Appendix F	Addition to the second-stage models

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Oxycodone hydrochloride*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Oxycodone hydroc...	12268908	21202	40	0914E-07734
Caelo	Oxycodone hydroc...	12034515	21204	40	0914F-07734
Caelo	Oxycodone hydroc...	12268914	22410	60	1416Q-07734
Caelo	Oxycodone hydroc...	12268913	22419	60	1416R-07734
Caelo	Oxycodone hydroc...	12268908	21202SI	60	0914E-07734

continued on the next page

continued from previous page

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Oxycodone hydroc...	12034515	21204SI	60	0914F-07734
Caelo	Oxycodone hydroc...	12268914	22410SI	60	1416Q-07734
Caelo	Oxycodone hydroc...	12268913	22419SI	60	1416R-07734

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 440 spectra of 8 reference samples from the substance/substance group *Oxycodone hydrochloride*. These samples are listed above in the [calibration samples](#) section.
- 107 809 spectra from a total of 1898 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 160 spectra of 4 reference samples from the substance/substance group *Oxycodone hydrochloride*.

Supplier	Substance	Batch	Sample ID	Spectra
Caelo	Oxycodone hydrochloride	15220205	23495	40
Caelo	Oxycodone hydrochloride	17270202	24661	40
Caelo	Oxycodone hydrochloride	15220205	23495SI	40
Caelo	Oxycodone hydrochloride	17270202	24661SI	40

- 242 806 spectra from a total of 3905 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 2 spectra from 1 *Apo-Ident* customers from 1 batches from the substance/substance group *Oxycodone hydrochloride*.

Supplier	Substance	Batch	Spectra
Caelo	Oxycodone hydrochloride	17270202	2

- 95 639 spectra from 646 *Apo-Ident* customers from a total of 30 194 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Oxycodone hydrochloride* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Oxycodone hydrochloride* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	440	0	107 809
Type B	0	160	0	242 806
Type C	0	2	0	95 639

The substance/substance group *Oxycodone hydrochloride* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 98.6364 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 96.2500 %)
Type C	100.0000 % (> 99.8794 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra (*Rule of Three* [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the *Mahalanobis distance* to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21202	21202	0.00	40.44
21204	21204	0.00	39.83
22410	22410	0.00	42.47
22419	22419	0.00	40.33

VALIDATION REPORT

IdentModul 2.0-2019-05

Validated substance/substance group **Phenobarbital**
Substance class Narcotic substances, solid
Report date 06/12/2019
Report number 20311-2019-12-06
Executing company HiperScan GmbH
Weißeritzstraße 3
01067 Dresden
Germany

Relevant substance names

Phenobarbital; Acidum phenylaethylbarbituricum

Special notes

When selecting the *Phenobarbital* substance/substance group, the following information is displayed to the user:

no information

Applicable documents

978-3-7692-6679-5 *European Pharmacopoeia 9th Edition, Basic Version 2017* [3]
Komm2.2.40 *Erfüllung von Ph. Eur. 2.2.40 durch Apo-Ident* [4]
AA004 *Erstellung und Validierung eines IdentModul-Updates*
Appendix F [Addition to the second-stage models](#)

Validation method

Validation is performed after every change to the chemometric model (also “database”) in three steps:

1. The chemometric model is calculated from the calibration spectra using a *PCA algorithm*. The calibration spectra originate from the calibration samples of all substances in this class.
2. In the *chemometric model* generated, the distances between all separable substances are checked for compliance with the specified safety distances.
3. All spectra are presented to the generated chemometric model for evaluation. In three runs, the reference spectra (*Type A*), spectra from independent samples (*Type B*), and spectra from the field (*Type C*) are presented successively. Here, no single *false positive* result is permitted.

Finally, a report is generated from the validation run results and is archived together with the parameters of the model generation in audit-compliant manner.

Calibration samples

Only spectra which have been recorded by *HiperScan GmbH* from traceable samples are used for the generation of the [chemometric models](#). The following samples have been obtained from the substance/substance group *Phenobarbital*:

Supplier	Substance	Batch	Sample ID	Spectra	Certificate
Caelo	Phenobarbital	12173202	21198	40	0914E-01040
Th. Geyer	Phenobarbital	67045	21200	40	0914F-01040
Fagron	Phenobarbital	15G07-B06-311718	22420	60	1416Q-01040
Th. Geyer	Phenobarbital	067056	23497	40	2017A-01040
Caelo	Phenobarbital	12173202	21198SI	60	0914E-01040
Th. Geyer	Phenobarbital	67045	21200SI	60	0914F-01040
Fagron	Phenobarbital	15G07-B06-311718	22420SI	60	1416Q-01040

Validation samples

A total of 446 856 spectra were provided for validation. The results were evaluated separately according to the following sample categories:

Type A All calibration spectra.

- 360 spectra of 7 reference samples from the substance/substance group *Phenobarbital*. These samples are listed above in the [calibration samples](#) section.
- 107 889 spectra from a total of 1898 batches from further 388 substances. These spectra were recorded by *HiperScan GmbH*. *Type A* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type A* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix A](#). The samples in [Appendix A](#) were added for model generation not with the intention to have the model identify these substances, but with the intention to make the model reject them despite their similarity to substances identifiable.

Type B Spectra from independent samples not included in the database structure. Measurement is carried out by *HiperScan GmbH*.

These samples are allowed to originate from the same batch as calibration samples (*Type A*) as long as the package differs. In the following table spectra recorded from the same package as the corresponding calibration samples (*Type A*) are highlighted by [†].

- 120 spectra of 3 reference samples from the substance/substance group *Phenobarbital*.

Supplier	Substance	Batch	Sample ID	Spectra
Fagron	Phenobarbital	16G20-B04-328344	23496	40
Fagron	Phenobarbital	16G20-B04-328344	23496SI	40
Th. Geyer	Phenobarbital	067056	23497SI	40

- 242 846 spectra from a total of 3905 batches from further 1077 substances. These spectra were recorded by *HiperScan GmbH*. *Type B* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type B* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix B](#).

Type C Spectra of independent samples not included in the database structure. *Apo-Ident* customers carry out the measurements. *HiperScan GmbH* generally does not verify the information provided by the customer regarding manufacturer and batch number.

- 7 spectra from 5 *Apo-Ident* customers from 4 batches from the substance/substance group *Phenobarbital*.

Supplier	Substance	Batch	Spectra
Th. Geyer	Phenobarbital	67045	2
Fagron	Phenobarbital	13G05-N09	1
BLS	Phenobarbital	67056	1
Th. Geyer	Phenobarbital	67056	2
Fagron	Phenobarbital	15G07-B06-311718	1

- 95 634 spectra from 646 *Apo-Ident* customers from a total of 30 191 batches from a further 976 substances. These spectra were recorded by *Apo-Ident* customers. *Type C* samples are documented within this entire report. If the substance can be identified with this [chemometric model](#), the samples are listed in the *Type C* section of each substance page. If the substance cannot be identified by this model, the samples are listed in [Appendix C](#).

Validation results

The validation runs checked whether the substance/substance group *Phenobarbital* can clearly be distinguished from all other substances using NIR spectroscopy with *Apo-Ident*. For this purpose, all relevant spectra of the various substances were compared with *Phenobarbital* and it was evaluated how many matches (positive) and rejections (negative) were correct or incorrect. The following table breaks down the numbers of correct and incorrect results according to the expected result (*positive/negative*) and the validation spectrum type (*A/B/C*).

	False positive	True positive	False negative	True negative
Type A	0	360	0	107 889
Type B	0	120	0	242 846
Type C	0	7	0	95 634

The substance/substance group *Phenobarbital* can be clearly distinguished from all other substances. In order to make these figures comparable, the weighted *true negative rate (specificity)* and the weighted *true positive rate (recognition rate)* are determined:

	Specificity	Recognition rate
Type A	100.0000 % (> 99.9908 %)	100.0000 % (> 98.3333 %)
Type B	100.0000 % (> 99.9935 %)	100.0000 % (> 95.0000 %)
Type C	100.0000 % (> 99.8792 %)	n/a (n/a)

In order to ensure that each substance is entered with the same weight irrespective of the number of spectra available, each spectrum is weighted with the reciprocal value of its number. (If several new spectra of a substance with a very distinctive spectrum are added, the values for *specificity* and *recognition rate* are safer, but are by no means idealised.)

In order to illustrate the impact of the number of spectra, a comparison is made in parentheses as to how the *specificity* or *recognition rate* would worsen if three additional incorrect results were to exist amongst the spectra ([Rule of Three](#) [10, 11]). The larger the number of spectra, the lower the deterioration if the three hypothetical *incorrect results* were added.

Where the number of spectra does not exceed 20, no *recognition rate* is provided.

Identity according to identity test using NIR spectroscopy

For each substance, a certificate of the correct identity is obtained from an independent GMP-certified test laboratory. If the identity of the sample can be provided using NIR on an independently tested reference sample, in the following table the [Mahalanobis distance](#) to the reference sample is specified as well as the *Mahalanobis distance* to the next non-identical substance:

Sample ID	Reference sample ID	Distance to reference sample	Distance to next foreign sample
21198	21198	0.00	49.27
21200	21200	0.00	50.03
22420	22420	0.00	53.01
23497	23497	0.00	51.69

Appendix A: Additional calibration samples (*Type A*)

Substance (spectra also enter the generation of the [chemometric model](#) which cannot be identified with this model. This serves the purpose that substances which the model is not supposed to recognise are safely rejected. The spectra of these samples are entered into the validation as *Type A*. They were recorded by *HiperScan GmbH*.

The samples originate from 2185 batches. From these, 101 298 spectra were recorded. The calibration spectra of substances which can be identified with the model are listed respectively in the section *calibration samples* for the individual substances and do not appear again elsewhere in this list.

Supplier	Substance	Batch	Spectra	Certificate
apomix	Non-ionic emulsi...	NS-8350010618	40	not necessary
apomix	Non-ionic emulsi...	NS-8350010216	40	not necessary
Audor	Alfatradiol	AP33-16-I-62P2985A	40	not necessary
Audor	Alfatradiol	AP33-18-I-68P3620A	40	not necessary
Audor	Benzydamine hydr...	AP300111801004B	40	not necessary
Audor	Colistin sulfate	APKB-CS-1602026C	40	not necessary
Audor	Colistin sulfate	APKB-CS-1603034C	40	not necessary
Audor	Colistin sulfate	APKB-CS-1602026C	40	not necessary
Audor	Cyproterone acet...	APX03-181002-m0C	40	not necessary
Audor	Cyproterone acet...	APCPA140201CP	40	not necessary
Audor	Hyaluronic acid ...	AP30101702H	40	1802037
Audor	Hyaluronic acid ...	AP30101404H	60	not necessary
Audor	Metoprolol tartr...	APMTP02218MT	40	not necessary
Audor	Naloxone hydroch...	AP1547I022-1N	40	not necessary
Audor	Naloxone hydroch...	AP13501078IIN	40	not necessary
Audor	Nystatin	AP4018717N	40	not necessary
Audor	Oxetacaine	AP010-160120	40	1809415
Audor	Prasterone	APZDHAC18004D	40	1906523
Audor	Prasterone	APZDHA16003-M	40	1805404
Audor	Prasterone	APZDHAC17003MD	40	not necessary
Audor	Tobramycin sulfate	AP15061047002T	40	1608462
Audor	Tobramycin sulfate	AP15061047002T	40	not necessary
Audor	Tranexamic acid	APX1803027MT	40	not necessary
Audor	Tranexamic acid	APX1508016MT	40	not necessary
Audor	Uridine	APIB160018211AU	40	not necessary
Audor	Vancomycin hydro...	APHAP1609006EV	40	AR-18-FG-000670-01
Auropure Life...	Amygdalin, bitte...	PA02099-17102131	40	from supplier
Bombastus	Hard fat in past...	308558	40	not necessary
Bombastus	Hard fat in past...	306363	34	not necessary
Bombastus	Non-ionic emulsi...	19001210	40	not necessary
Bombastus	Non-ionic emulsi...	19001210	40	not necessary
Caelo	5-aminolevulinic...	16357904	40	1707616
Caelo	5-aminolevulinic...	190219	40	20190312*
Caelo	5-aminolevulinic...	16357904	40	not necessary
Caelo	Acetylcysteine	12299020	40	not necessary
Caelo	Acetylcysteine	161398	60	20160520*
Caelo	Acetylsalicylic ...	13297509	40	not necessary
Caelo	Acetylsalicylic ...	161316	40	20160825*
Caelo	Alfatradiol	13201601	120	not necessary
Caelo	Alfatradiol	191599	40	20190812*
Caelo	Alfatradiol	13201601	60	not necessary
Caelo	Allantoin	172346	40	20171019*
Caelo	Allantoin	162289	40	20160815*
Caelo	Allantoin	172346	40	20171019*
Caelo	Ambroxol hydroch...	171582	40	20170607*
Caelo	Ambroxol hydroch...	161996	40	20160810*
Caelo	Amifampridine	17251208	40	not necessary
Caelo	Ascorbic acid	171358	40	20170524*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Ascorbic acid	171358	40	20170524*
Caelo	Ascorbic acid	12350011	40	not necessary
Caelo	Ascorbic acid	163452	40	20170125*
Caelo	Atropine sulfate	15301601	60	not necessary
Caelo	Atropine sulfate	14336806	60	1510209
Caelo	Atropine sulfate	14336806	60	1510209
Caelo	Atropine sulfate	11180602	30	1405538
Caelo	Atropine sulfate	13173808	40	1403483
Caelo	Atropine sulfate	15301601	60	not necessary
Caelo	Atropine sulfate	13357405	50	1504594
Caelo	Benzocaine	13213509	60	1706017
Caelo	Benzocaine	12341904	40	1411055
Caelo	Benzoic acid	15084614	40	not necessary
Caelo	Benzoic acid	13391310	60	not necessary
Caelo	Benzoic acid	15084614	40	not necessary
Caelo	Benzoic acid	15084614	40	1707644
Caelo	Benzoic acid	15084611	40	not necessary
Caelo	Betaine monohydr...	12336602	40	1403327
Caelo	Betamethasone di...	13372501	45	AR-15-FG-005428-01
Caelo	Betamethasone di...	13372501	40	AR-15-FG-005428-01
Caelo	Betamethasone di...	13236803	46	not necessary
Caelo	Betamethasone di...	13236803	50	not necessary
Caelo	Betamethasone di...	13236803	40	not necessary
Caelo	Betamethasone va...	13225507	54	1706018
Caelo	Betamethasone va...	13091213	40	not necessary
Caelo	Bifonazole	15152701	40	not necessary
Caelo	Bifonazole	15152702	40	not necessary
Caelo	Bifonazole	15152702	40	not necessary
Caelo	Bifonazole	151527	40	20150703*
Caelo	Bilberry fruit, ...	10265602	40	not necessary
Caelo	Biotin	16164313	40	not necessary
Caelo	Biotin	163432	40	20170119*
Caelo	Biotin	16164313	40	not necessary
Caelo	Biotin	163432	40	20170119*
Caelo	Bismuth gallath,...	13179105	40	1412254
Caelo	Bismuth gallath,...	170041	40	20170214*
Caelo	Bismuth gallath,...	17004102	40	not necessary
Caelo	Bismuth gallath,...	13179105	60	not necessary
Caelo	Black mustard se...	12350204	80	1403711
Caelo	Boric acid	14133119	60	not necessary
Caelo	Boric acid	170613	40	20170316*
Caelo	Boric acid	15425112	40	not necessary
Caelo	Budesonide, micr...	15344004	40	not necessary
Caelo	Budesonide, micr...	11351704	40	AR-15-FG-005286-01
Caelo	Budesonide, micr...	161431	25	20160606*
Caelo	Budesonide, micr...	15344004	23	not necessary
Caelo	Budesonide, micr...	161431	40	20160606*
Caelo	Butylated hydrox...	15172103	60	not necessary
Caelo	Butylated hydrox...	15172103	60	not necessary
Caelo	Butylated hydrox...	13354604	40	AR-15-FG-007846-02
Caelo	Butylated hydrox...	15172103	60	not necessary
Caelo	Butylated hydrox...	13354602	40	not necessary
Caelo	Butylated hydrox...	190749	40	20190405*
Caelo	Caffeine	13233105	60	not necessary
Caelo	Caffeine	13233104	60	not necessary
Caelo	Caffeine	12317911	40	not necessary
Caelo	Caffeine	12317911	40	not necessary
Caelo	Calcium ascorbate	16311101	40	not necessary
Caelo	Calcium ascorbate	163111	40	20161207*
Caelo	Calcium ascorbate	13105624	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Calcium citrate ...	13391203	40	not necessary
Caelo	Calcium citrate ...	13391203	40	not necessary
Caelo	Calcium citrate ...	171768	40	20170713*
Caelo	Calcium gluconate	12335829	40	not necessary
Caelo	Calcium gluconate	12236502	40	1401575
Caelo	Calcium glycerop...	14341007	60	AR-16-FG-005932-01
Caelo	Calcium glycerop...	13007912	40	not necessary
Caelo	Calcium glycerop...	191007	40	20190502*
Caelo	Calcium lactate ...	161313	60	20160513*
Caelo	Calcium lactate ...	15101401	60	not necessary
Caelo	Calcium lactate ...	15101405	60	not necessary
Caelo	Calcium pantothe...	14208012	60	not necessary
Caelo	Calcium pantothe...	173352	40	20180202*
Caelo	Camphor, racemic	15316801	60	not necessary
Caelo	Camphor, racemic	170043	40	20170220*
Caelo	Captopril	12339610	40	1403227
Caelo	Captopril	190118	40	20190312*
Caelo	Captopril	15189502	60	AR-15-FG-014855-01
Caelo	Carbomers 35000 ...	14057902	40	not necessary
Caelo	Carbomers 35000 ...	13135601	40	1403767
Caelo	Carbomers 35000 ...	13186207	40	1403768
Caelo	Carbomers 35000 ...	14057903	60	not necessary
Caelo	Carbomers 35000 ...	162530	40	20161011*
Caelo	Carmellose sodiu...	13299310	60	not necessary
Caelo	Carmellose sodiu...	14029211	40	not necessary
Caelo	Carmellose sodiu...	16171406	40	AR-17-FG-015964-01
Caelo	Carmellose sodiu...	13104807	40	1401186
Caelo	Carmellose sodiu...	180486	40	20180228*
Caelo	Carmellose sodiu...	12354404	40	1401185
Caelo	Carmellose sodiu...	161714	40	20160713*
Caelo	Cetomacrogol emu...	14055001	40	not necessary
Caelo	Cetomacrogol emu...	15277504	40	not necessary
Caelo	Cetostearyl alco...	14054504	60	not necessary
Caelo	Cetostearyl alco...	17049205	40	not necessary
Caelo	Cetostearyl alco...	14054504	40	not necessary
Caelo	Cetyl alcohol	161063	40	20160420*
Caelo	Cetyl alcohol	161063	40	20160420*
Caelo	Cetyl alcohol	161063	40	20160420*
Caelo	Cetyl alcohol	161063	40	20160420*
Caelo	Cetyl alcohol	161063	40	20160420*
Caelo	Cetyl alcohol	13286502	60	not necessary
Caelo	Cetyl alcohol	15323301	60	AR-16-FG-008835-01
Caelo	Cetyl alcohol	161063	40	20160420*
Caelo	Cetyl alcohol	13096901	40	1405179
Caelo	Cetyl alcohol	13286502	40	not necessary
Caelo	Cetyl palmitate 15	160538	60	20160323*
Caelo	Cetyl palmitate 15	191390	40	20190702*
Caelo	Cetyl palmitate 15	172715	40	20171124*
Caelo	Cetyl palmitate 15	14351606	60	not necessary
Caelo	Cetyl palmitate 15	171640	40	20170719*
Caelo	Cetyl palmitate 15	15301703	60	AR-16-FG-007435-01
Caelo	Cetyl palmitate 15	17043503	40	not necessary
Caelo	Cetyl palmitate 15	15301703	60	not necessary
Caelo	Cetyl palmitate 15	12331916	40	1402429
Caelo	Chicory root, cut	11374504	40	not necessary
Caelo	Chloral hydrate	13213311	40	not necessary
Caelo	Chloral hydrate	15448318	40	1704462
Caelo	Chloral hydrate	16355802	40	not necessary
Caelo	Chloramphenicol	13080924	60	1509287
Caelo	Chloramphenicol	11370716	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Chloramphenicol	12248507	40	not necessary
Caelo	Chloramphenicol	14346908	60	not necessary
Caelo	Chlorhexidine di...	15438909	40	20160309*
Caelo	Chlorhexidine di...	154389	60	20160309*
Caelo	Chlorhexidine di...	154389	40	20160309*
Caelo	Chlorhexidine di...	154389	40	20160309*
Caelo	Chlorhexidine di...	15438913	40	20160309*
Caelo	Chlorhexidine di...	12142726	40	1402685
Caelo	Chlorhexidine di...	15438914	40	1704392
Caelo	Chlorhexidine di...	15438914	40	not necessary
Caelo	Chlorhexidine di...	154389	40	20160309*
Caelo	Chlorhexidine di...	12142727	40	not necessary
Caelo	Chlorocresol	163378	40	20170124*
Caelo	Chlorocresol	13275106	60	not necessary
Caelo	Chlorocresol	182523	40	20180919*
Caelo	Chlorocresol	13275106	60	not necessary
Caelo	Chlorocresol	190837	40	20190508*
Caelo	Cholesterol	16044415	40	1804387
Caelo	Cholesterol	160444	60	20160222*
Caelo	Choline chloride	162092	40	20160801*
Caelo	Choline chloride	14273819	40	not necessary
Caelo	Choline chloride	14273819	40	not necessary
Caelo	Choline chloride	16209209	40	not necessary
Caelo	Ciclopirox olamine	14147406	60	AR-15-FG-014857-01
Caelo	Ciclopirox olamine	13343209	40	AR-15-FG-016214-01
Caelo	Ciclopirox olamine	14378809	60	not necessary
Caelo	Ciclopirox olamine	15364408	40	not necessary
Caelo	Ciclopirox olamine	15364408	40	not necessary
Caelo	Ciclopirox olamine	14378809	60	not necessary
Caelo	Ciclopirox olamine	14147406	60	AR-15-FG-014857-01
Caelo	Ciclosporin	173327	35	20180111*
Caelo	Ciclosporin	13148604	60	not necessary
Caelo	Ciclosporin	13148604	60	1412541
Caelo	Citric acid	11380204	40	not necessary
Caelo	Citric acid	13396502	60	not necessary
Caelo	Citric acid	160386	60	20160315*
Caelo	Citric acid	13243213	40	not necessary
Caelo	Citric acid	12043112	40	not necessary
Caelo	Citric acid mono...	15185301	60	not necessary
Caelo	Citric acid mono...	17085602	40	not necessary
Caelo	Citric acid mono...	181305	40	20180606*
Caelo	Clindamycin hydr...	15392817	40	not necessary
Caelo	Clindamycin hydr...	15392809	40	not necessary
Caelo	Clindamycin-2-di...	163265	40	20161214*
Caelo	Clindamycin-2-di...	13166605	40	1403485
Caelo	Clindamycin-2-di...	13166605	40	not necessary
Caelo	Clindamycin-2-di...	16326504	40	AR-17-FG-022557-01
Caelo	Clindamycin-2-di...	15435801	60	AR-16-FG-009129-01
Caelo	Clindamycin-2-di...	191224	40	20190523*
Caelo	Clioquinol	11304120	40	not necessary
Caelo	Clioquinol	163072	40	20161202*
Caelo	Clobetasol propi...	13134705	80	AR-15-FG-005328-01
Caelo	Clobetasol propi...	15439504	40	not necessary
Caelo	Clobetasol propi...	14338001	59	AR-16-FG-008581-01
Caelo	Clobetasol propi...	14338001	59	1607072
Caelo	Clobetasol propi...	15439504	37	not necessary
Caelo	Clobetasol propi...	14338001	60	AR-16-FG-008581-01
Caelo	Clobetasol propi...	14338001	45	1607072
Caelo	Clobetasol propi...	13269505	60	AR-15-FG-000622-02
Caelo	Clotrimazole	13224204	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Clotrimazole	13224220	40	not necessary
Caelo	Clotrimazole	13224220	40	not necessary
Caelo	Clotrimazole	13184113	40	not necessary
Caelo	Clotrimazole	153956	40	20151215*
Caelo	Cocoa butter	17055907	33	not necessary
Caelo	Cocoa butter	190360	40	20190227*
Caelo	Cocoa butter	14348906	60	not necessary
Caelo	Cocoa butter	163165	40	20161116*
Caelo	Cocoa butter	13292013	40	not necessary
Caelo	Colophony	191664	40	20190701*
Caelo	Colophony	12301410	60	AR-15-FG-016314-01
Caelo	Colophony	161332	40	20160510*
Caelo	Colophony	192081	40	20190731*
Caelo	Colophony	15186804	60	AR-16-FG-005931-01
Caelo	Corn starch	13154914	60	not necessary
Caelo	Corn starch	182748	40	20181031*
Caelo	Corn starch	161385	60	20160621*
Caelo	Cyanocobalamin	13342803	40	AR-15-FG-005429-01
Caelo	Cyanocobalamin	160154	60	20160309*
Caelo	Cyanocobalamin	13342807	60	AR-15-FG-014917-01
Caelo	Cyanocobalamin	160154	60	20160309*
Caelo	Cyanocobalamin	13342807	30	AR-15-FG-014917-01
Caelo	Dapsone	160735	60	20160405*
Caelo	Dapsone	191135	40	20190726*
Caelo	Dapsone	170311	39	20170221*
Caelo	Dequalinium chlo...	15402803	40	1707614
Caelo	Dequalinium chlo...	15218406	40	not necessary
Caelo	Dequalinium chlo...	154028	40	20160219*
Caelo	Dequalinium chlo...	154028	40	20160219*
Caelo	Dequalinium chlo...	12361908	40	1412076
Caelo	Dequalinium chlo...	12361906	60	not necessary
Caelo	Dexamethasone	190221	40	20190214*
Caelo	Dexamethasone	15054627	40	not necessary
Caelo	Dexamethasone	13084315	60	AR-15-FG-006810-01
Caelo	Dexamethasone	190221	40	20190214*
Caelo	Dexamethasone ac...	14200207	60	not necessary
Caelo	Dexamethasone ac...	14200207	30	not necessary
Caelo	Dexamethasone ac...	13260706	60	1502537
Caelo	Dexamethasone ac...	13260703	40	1412430
Caelo	Diclofenac sodium	12214219	20	not necessary
Caelo	Diclofenac sodium	13214311	60	not necessary
Caelo	Diclofenac sodium	163492	40	20170112*
Caelo	Diclofenac sodium	13214311	60	not necessary
Caelo	Diltiazem hydroc...	13346914	60	not necessary
Caelo	Diltiazem hydroc...	13346914	60	not necessary
Caelo	Diltiazem hydroc...	13346914	40	not necessary
Caelo	Diltiazem hydroc...	161908	40	20160715*
Caelo	Dimethyl fumarate	13269801	60	1506634
Caelo	Diphenhydramine ...	16004901	40	not necessary
Caelo	Diphenhydramine ...	12329816	60	not necessary
Caelo	Diphenhydramine ...	180082	40	20180126*
Caelo	Diphenhydramine ...	12329816	40	not necessary
Caelo	Dithranol	18190604	40	1904372
Caelo	Dithranol	190868	40	20190426*
Caelo	Dithranol	190550	40	20190426*
Caelo	Dithranol	18115302	40	1811302
Caelo	Doxylamine hydro...	183312	40	20190219*
Caelo	Doxylamine hydro...	15203906	40	not necessary
Caelo	Doxylamine hydro...	15203906	40	not necessary
Caelo	Emulsifying ceto...	13233005	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Emulsifying ceto...	153857	60	20151202*
Caelo	Emulsifying ceto...	153856	60	20160210*
Caelo	Emulsifying ceto...	153857	60	20151202*
Caelo	Emulsifying ceto...	14322505	60	not necessary
Caelo	Emulsifying ceto...	153856	60	20160210*
Caelo	Emulsifying ceto...	13233005	30	1506057
Caelo	Ephedrine hydroc...	13323305	40	not necessary
Caelo	Ephedrine hydroc...	14378502	60	not necessary
Caelo	Ephedrine hydroc...	161197	60	20160420*
Caelo	Ephedrine hydroc...	12017217	40	1511635
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	60	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	60	1603180
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	60	1603180
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	60	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	12337804	40	not necessary
Caelo	Epinephrine hydr...	12120202	40	1405614
Caelo	Erythromycin	160792	60	20160408*
Caelo	Erythromycin	13175601	40	not necessary
Caelo	Estradiol benzoate	181444	40	20180614*
Caelo	Estradiol benzoate	15078608	40	not necessary
Caelo	Estradiol hemihy...	13285802	60	not necessary
Caelo	Estradiol hemihy...	13285802	60	1510622
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	60	1410504
Caelo	Estradiol hemihy...	12268016	40	1410502
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	163277	40	20161221*
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	13264001	120	1510624
Caelo	Estradiol valerate	13264001	60	1410503
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	13264001	60	1410503
Caelo	Estradiol valerate	15194802	60	not necessary
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	163277	40	20161221*
Caelo	Estradiol valerate	13264001	20	not necessary
Caelo	Estradiol valerate	15194802	30	not necessary
Caelo	Estradiol valerate	13264001	45	1510624
Caelo	Estradiol valerate	13264001	45	1410503
Caelo	Estriol	13256001	60	1511155
Caelo	Estriol	13256001	55	1511155
Caelo	Estriol	13256001	60	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	60	not necessary
Caelo	Estriol	13256001	20	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	12298706	40	not necessary
Caelo	Ethacridine lact...	12149507	40	not necessary
Caelo	Ethacridine lact...	16273403	40	not necessary
Caelo	Fampridine	180219	40	20180214*
Caelo	Fampridine	190675	40	20190506*
Caelo	Fampridine	180219	40	20180214*
Caelo	Fampridine	16276001	40	not necessary
Caelo	Fampridine	16276001	40	not necessary
Caelo	Fampridine	16276001	40	not necessary
Caelo	Fampridine	16276001	40	not necessary
Caelo	Fluocinolonaceto...	12035003	48	not necessary
Caelo	Fluocinolonaceto...	162281	40	20160907*
Caelo	Fluocinolonaceto...	12035003	60	not necessary
Caelo	Fluocinolonaceto...	162281	30	20160907*
Caelo	Fructose	160667	60	20160425*
Caelo	Fructose	15084311	40	not necessary
Caelo	Gelatine, white	15320504	40	AR-16-FG-014461-01
Caelo	Gelatine, white	15320504	45	AR-16-FG-011564-01
Caelo	Gelatine, white	11334025	40	not necessary
Caelo	Gentamicin sulfate	13205003	60	not necessary
Caelo	Gentamicin sulfate	14353504	40	not necessary
Caelo	Ginger, peeled	11325715	40	not necessary
Caelo	Glucose	14033901	40	not necessary
Caelo	Glucose	14142310	45	not necessary
Caelo	Glucose	16050501	40	not necessary
Caelo	Glucose	12294514	40	1402475
Caelo	Glucose	190380	40	20190411*
Caelo	Glucose monohydr...	173050	40	20171110*
Caelo	Glucose monohydr...	171368	40	20170601*
Caelo	Glycerol monoste...	14099903	60	AR-15-FG-000732-02
Caelo	Glycerol monoste...	12105705	80	not necessary
Caelo	Glycerol monoste...	154450	60	20160225*
Caelo	Glycerol monoste...	12105720	40	not necessary
Caelo	Glycerol monoste...	15286302	60	not necessary
Caelo	Glycerol monoste...	182462	40	20181008*
Caelo	Glycerol monoste...	13232604	60	1504599
Caelo	Glycerol monoste...	190547	40	20190423*
Caelo	Glycerol monoste...	180683	40	20180430*
Caelo	Glycerol monoste...	13109304	40	1407633
Caelo	Glycerol monoste...	16233804	40	10052
Caelo	Glycerol monoste...	15286302	60	AR-16-FG-009226-01
Caelo	Glycerol monoste...	12210401	40	not necessary
Caelo	Glycerol monoste...	12210401	80	1407214
Caelo	Glycerol monoste...	161646	40	20160727*
Caelo	Glycerol monoste...	16030301	40	not necessary
Caelo	Glycerol monoste...	14240607	60	not necessary
Caelo	Glycerol monoste...	14240607	60	1412120
Caelo	Glycine	170265	40	20170315*
Caelo	Glycine	13157303	60	not necessary
Caelo	Glycine	15262518	40	not necessary
Caelo	Glycine	14147006	60	AR-16-FG-012263-01
Caelo	Glycine	173493	40	20180117*
Caelo	Guaifenesin	172496	40	20171005*
Caelo	Guaifenesin	163333	40	20161221*
Caelo	Hard fat in past...	12158122	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Hard fat in past...	160278	60	20160407*
Caelo	Hay flower	12317104	40	not necessary
Caelo	Hydrochlorothiaz...	182649	40	20181127*
Caelo	Hydrochlorothiaz...	13008505	60	AR-15-FG-014461-01
Caelo	Hydrocortisone	180915	40	20180412*
Caelo	Hydrocortisone	13341502	60	not necessary
Caelo	Hydrocortisone	13341502	60	not necessary
Caelo	Hydrocortisone a...	182780	40	20181210*
Caelo	Hydrocortisone a...	13291708	60	not necessary
Caelo	Hydroquinone	14281501	60	AR-15-FG-008058-01
Caelo	Hydroquinone	14026904	60	not necessary
Caelo	Hydroquinone	14280110	60	not necessary
Caelo	Hydroquinone	12288505	40	not necessary
Caelo	Hydroquinone	182306	40	20181008*
Caelo	Hydroquinone	15057001	60	not necessary
Caelo	Hydroquinone	15057004	60	not necessary
Caelo	Hydroquinone	15057001	60	not necessary
Caelo	Hydroxyethylcell...	12249204	40	not necessary
Caelo	Hydroxyethylcell...	16278705	40	not necessary
Caelo	Hydroxyethylcell...	161006	60	20160426*
Caelo	Hydroxyethylcell...	16278705	40	not necessary
Caelo	Hydroxyethylcell...	13378506	40	1507205
Caelo	Hydroxyethylcell...	190769	40	20190425*
Caelo	Hydroxyethylcell...	171151	40	20170523*
Caelo	Hydroxyethylcell...	16252008	40	not necessary
Caelo	Hydroxyethylcell...	14352504	60	AR-15-FG-009318-01
Caelo	Hydroxyethylcell...	03936475	20	not necessary
Caelo	Hydroxypropylcel...	181146	40	20180618*
Caelo	Hydroxypropylcel...	14189804	60	AR-15-FG-014922-01
Caelo	Hydroxypropylcel...	12093107	40	1402049
Caelo	Hydroxypropylcel...	15345001	60	not necessary
Caelo	Hyoscine hydrobr...	15054501	28	1511156
Caelo	Hyoscine hydrobr...	14278801	60	1504154
Caelo	Hyoscine hydrobr...	15054501	60	1511156
Caelo	Hypromellose 4000	14338801	45	not necessary
Caelo	Hypromellose 4000	13155006	40	1402056
Caelo	Hypromellose 4000	13155006	40	not necessary
Caelo	Ibuprofen	172271	40	20170831*
Caelo	Ibuprofen	11352219	40	not necessary
Caelo	Iodoform	16012302	40	not necessary
Caelo	Iodoform	191132	40	20190524*
Caelo	Iodoform	12283206	40	not necessary
Caelo	Iodoform	16012302	40	AR-17-FG-012459-01
Caelo	Isoxuprine hydro...	16177001	40	not necessary
Caelo	Isoxuprine hydro...	161770	40	20160818*
Caelo	Isoxuprine hydro...	181078	40	20180719*
Caelo	Ispaghula husk	13238005	40	not necessary
Caelo	Juniper, pounded	13367302	40	not necessary
Caelo	Kelp, cut	15418801	60	AR-17-FG-007505-01
Caelo	Kelp, cut	13230805	40	AR-15-FG-007678-01
Caelo	Ketoconazole	13179612	60	not necessary
Caelo	Ketoconazole	12171007	30	1401554
Caelo	Ketoconazole	13179612	60	not necessary
Caelo	Ketoconazole	153477	40	20151113*
Caelo	Lactose monohydr...	14025601	40	not necessary
Caelo	Lactose monohydr...	14025601	40	not necessary
Caelo	Lactose monohydr...	12139910	30	not necessary
Caelo	Lactose monohydr...	14025601	40	not necessary
Caelo	Lactose monohydr...	14025601	40	not necessary
Caelo	Lactose monohydr...	14025601	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Lactose monohydr...	160277	60	20160321*
Caelo	Lidocaine	180444	40	20180419*
Caelo	Lidocaine	13352411	60	not necessary
Caelo	Lidocaine hydroc...	172923	40	20171122*
Caelo	Lidocaine hydroc...	12135511	40	AR-15-FG-008867-01
Caelo	Liquorice root, ...	11351911	40	not necessary
Caelo	L-phenylalanine	14088304	60	1501656
Caelo	Macrogols 1500	17120001	40	not necessary
Caelo	Macrogols 1500	12324607	40	1403323
Caelo	Macrogols 1500	171200	40	20170519*
Caelo	Macrogols 4000	12361502	40	not necessary
Caelo	Macrogols 4000	13356616	40	not necessary
Caelo	Macrogols 4000	12074116	40	1401164
Caelo	Macrogols 4000	12318019	40	1403322
Caelo	Magnesium hydrog...	15415014	40	not necessary
Caelo	Magnesium hydrog...	13389102	60	not necessary
Caelo	Magnesium orotat...	16222701	40	not necessary
Caelo	Magnesium orotat...	162227	40	20160809*
Caelo	Magnesium stearate	13095536	40	not necessary
Caelo	Magnesium stearate	14136929	60	AR-15-FG-012028-01
Caelo	Mannitol	15197305	40	not necessary
Caelo	Mannitol	12305723	40	1402329
Caelo	Marshmallow root...	13214804	40	not necessary
Caelo	Matricaria flower	12337605	40	not necessary
Caelo	Melatonin	12321411	40	not necessary
Caelo	Melatonin	12321404	40	not necessary
Caelo	Melatonin	12321415	80	not necessary
Caelo	Menthol	13377404	60	not necessary
Caelo	Menthol	181522	40	20180625*
Caelo	Metamizol sodium...	15362816	40	not necessary
Caelo	Metamizol sodium...	161610	40	20160714*
Caelo	Methacholine chl...	16229304	40	not necessary
Caelo	Methacholine chl...	16229301	40	not necessary
Caelo	Methacholine chl...	13260303	40	not necessary
Caelo	Methacholine chl...	162293	40	20160815*
Caelo	Methacholine chl...	13260303	60	not necessary
Caelo	Methionine	16274007	40	not necessary
Caelo	Methionine	15404102	40	AR-17-FG-000285-02
Caelo	Methoxsalen	170781	40	20170427*
Caelo	Methoxsalen	16103003	40	not necessary
Caelo	Methoxsalen	170781	40	20170427*
Caelo	Methyl-4-hydroxy...	11090034	40	not necessary
Caelo	Methyl-4-hydroxy...	181769	40	20180824*
Caelo	Methylcellulose	13130115	60	not necessary
Caelo	Methylcellulose	160602	60	20160323*
Caelo	Metoclopramide h...	14193102	60	AR-15-FG-002889-01
Caelo	Metoclopramide h...	15421101	40	not necessary
Caelo	Metoclopramide h...	14193102	60	not necessary
Caelo	Metronidazole	13095007	40	not necessary
Caelo	Metronidazole	13095007	40	not necessary
Caelo	Metronidazole	160892	60	20160413*
Caelo	Metronidazole	12110410	40	not necessary
Caelo	Miconazole nitrate	13185601	80	not necessary
Caelo	Miconazole nitrate	161447	40	20160617*
Caelo	Miconazole nitrate	13352802	40	not necessary
Caelo	Miconazole nitrate	13185601	40	not necessary
Caelo	Miconazole nitrate	13352802	39	1412250
Caelo	Minoxidil	13039501	40	not necessary
Caelo	Minoxidil	14093701	40	not necessary
Caelo	Minoxidil	14093701	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Mometasone furoa...	163532	40	20170104*
Caelo	Mometasone furoa...	160286	60	20160316*
Caelo	Mupirocin calcium	14276602	60	not necessary
Caelo	Mupirocin calcium	161184	60	20160504*
Caelo	Naloxone hydroch...	163577	40	20170306*
Caelo	Naloxone hydroch...	170441	40	20170306*
Caelo	Naphazoline hydr...	181193	40	20180504*
Caelo	Naphazoline hydr...	12375211	60	not necessary
Caelo	Naphazoline hydr...	12375211	60	not necessary
Caelo	Naphazoline hydr...	12375211	40	not necessary
Caelo	Neomycin sulfate	161650	30	20160629*
Caelo	Neomycin sulfate	13062923	60	AR-15-FG-016581-01
Caelo	Neomycin sulfate	13062923	45	not necessary
Caelo	Neostigmine brom...	14256802	60	AR-15-FG-014460-01
Caelo	Neostigmine brom...	12172304	40	1403163
Caelo	Neostigmine brom...	10144105	40	not necessary
Caelo	Neostigmine brom...	15228101	60	not necessary
Caelo	Neostigmine brom...	14256802	60	AR-15-FG-014460-01
Caelo	Neostigmine brom...	15228101	60	not necessary
Caelo	Nicotinamide	163048	40	20161123*
Caelo	Nicotinamide	13097506	40	not necessary
Caelo	Nicotinamide	163048	40	20161123*
Caelo	Non-ionic emulsi...	13413001	40	not necessary
Caelo	Non-ionic emulsi...	13413001	40	1405257
Caelo	Non-ionic emulsi...	182021	40	20180725*
Caelo	Norethisterone a...	15417004	40	1711489
Caelo	Norethisterone a...	182252	40	20181031*
Caelo	Nystatin	17025110	40	not necessary
Caelo	Nystatin	16272214	40	not necessary
Caelo	Octenidine dihyd...	15275502	60	not necessary
Caelo	Octenidine dihyd...	13071407	40	1407345
Caelo	Octenidine dihyd...	13166806	40	not necessary
Caelo	Octenidine dihyd...	15356601	60	1604116
Caelo	Octenidine dihyd...	18223413	40	not necessary
Caelo	Octenidine dihyd...	163355	40	20170403*
Caelo	Octenidine dihyd...	15356601	60	1604116
Caelo	Octenidine dihyd...	170423	40	20170214*
Caelo	Octenidine dihyd...	12115202	40	not necessary
Caelo	Octenidine dihyd...	18223413	40	not necessary
Caelo	Octenidine dihyd...	170423	40	20170214*
Caelo	Orris root	11026705	40	not necessary
Caelo	Oxalic acid	171140	40	20170426*
Caelo	Oxalic acid	171140	40	20170426*
Caelo	Oxalic acid	13353901	40	not necessary
Caelo	Oxalic acid	18202403	40	1905160
Caelo	Oxalic acid	13353918	60	not necessary
Caelo	Oxalic acid	12212410	40	1403223
Caelo	Oxetacaine	191134	40	20190523*
Caelo	Oxetacaine	172498	40	20170925*
Caelo	Oxytetracycline ...	16077002	40	AR-16-FG-011460-01
Caelo	Oxytetracycline ...	12328103	40	not necessary
Caelo	Oxytetracycline ...	12130001	60	not necessary
Caelo	Pancreas powder	16139708	40	not necessary
Caelo	Pancreas powder	16139708	40	not necessary
Caelo	Pancreas powder	161397	60	20160527*
Caelo	Pancreas powder	13118902	40	1402190
Caelo	Pancreas powder	13395203	60	not necessary
Caelo	Papaverine hydro...	123161	40	20121127*
Caelo	Papaverine hydro...	13324310	60	not necessary
Caelo	Paracetamol	15075707	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Paracetamol	15075706	40	not necessary
Caelo	Paraffin, hard	13390413	60	not necessary
Caelo	Paraffin, hard	15232004	60	1702507
Caelo	Paraffin, microc...	154128	60	20160112*
Caelo	Paraffin, microc...	14307106	60	not necessary
Caelo	Paraffin, microc...	14307104	60	AR-15-FG-005277-01
Caelo	Patent blue V	13351807	40	1702362
Caelo	Patent blue V	12364206	40	not necessary
Caelo	Patent blue V	12364206	40	1508043
Caelo	Patent blue V	190746	40	20190425*
Caelo	Patent blue V	13351803	60	not necessary
Caelo	Phenazone	15415107	40	not necessary
Caelo	Phenazone	181347	40	20180629*
Caelo	Phenazone	15415101	40	not necessary
Caelo	Pheniramine male...	12200502	40	1401559
Caelo	Pheniramine male...	15172702	40	1908273
Caelo	Phenylbutazon	161121	60	20160331*
Caelo	Phenylbutazon	14057202	40	not necessary
Caelo	Pilocarpine hydr...	12217319	60	not necessary
Caelo	Pilocarpine hydr...	171595	40	20170606*
Caelo	Pilocarpine hydr...	17124208	40	not necessary
Caelo	Pilocarpine hydr...	12217319	60	1501660
Caelo	Poloxamers 407	16242605	40	not necessary
Caelo	Poloxamers 407	14196101	60	AR-14-FG-013090-01
Caelo	Polymyxin B sulf...	153501	60	20160318*
Caelo	Polymyxin B sulf...	181154	40	20180507*
Caelo	Potassium citrate	11320212	40	not necessary
Caelo	Potassium citrate	15123103	60	not necessary
Caelo	Potassium citrate	13317618	60	not necessary
Caelo	Potassium citrate	161130	60	20160602*
Caelo	Potassium hydrog...	170656	40	20170511*
Caelo	Potassium hydrog...	15354609	40	not necessary
Caelo	Potassium hydrog...	11345805	40	not necessary
Caelo	Potassium hydrog...	182051	40	20180903*
Caelo	Potassium sodium...	12030207	40	not necessary
Caelo	Potassium sodium...	14376001	60	not necessary
Caelo	Potassium sodium...	13216605	40	1404530
Caelo	Potassium sorbate	163266	40	20170105*
Caelo	Potassium sorbate	17170008	30	not necessary
Caelo	Potassium sorbate	13224604	60	not necessary
Caelo	Potassium tartra...	13331104	60	not necessary
Caelo	Potassium tartra...	15167701	60	not necessary
Caelo	Potassium tartra...	160289	60	20160405*
Caelo	Potato starch	161237	60	20160621*
Caelo	Potato starch	12360213	40	1403494
Caelo	Potato starch	14033815	60	not necessary
Caelo	Povidone 25, sol...	14142103	40	AR-15-FG-000596-01
Caelo	Povidone 25, sol...	16307808	40	AR-17-FG-015299-02
Caelo	Povidone 25, sol...	191752	40	20190717*
Caelo	Povidone 25, sol...	191276	40	20190517*
Caelo	Povidone 25, sol...	180276	40	20180202*
Caelo	Povidone 25, sol...	12032714	40	1402437
Caelo	Povidone iodine	15234101	40	AR-17-FG-015433-02
Caelo	Povidone iodine	13337205	60	not necessary
Caelo	Povidone iodine	182385	40	20180917*
Caelo	Povidone iodine	15234109	40	not necessary
Caelo	Prednicarbate, m...	13372404	33	not necessary
Caelo	Prednicarbate, m...	13119402	30	not necessary
Caelo	Prednicarbate, m...	13119402	60	not necessary
Caelo	Prednicarbate, m...	13119402	60	1501576

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Prednicarbate, m...	13372404	49	not necessary
Caelo	Prednicarbate, m...	13372404	50	AR-15-FG-007941-02
Caelo	Prednisolon dihy...	13304704	60	1411052
Caelo	Prednisolon dihy...	13304703	40	not necessary
Caelo	Prednisolon dihy...	16174715	40	not necessary
Caelo	Prednisolone ace...	12091814	80	not necessary
Caelo	Prednisolone ace...	13262605	47	not necessary
Caelo	Prednisolone ace...	13262605	60	not necessary
Caelo	Prednisolone ace...	13262605	60	not necessary
Caelo	Prednisolone ace...	12324502	40	not necessary
Caelo	Prednisolone, mi...	13343105	60	1411232
Caelo	Prednisolone, mi...	161343	60	20160601*
Caelo	Prednisolone, mi...	15400002	40	not necessary
Caelo	Prednisolone, mi...	15400002	40	not necessary
Caelo	Prednisolone, mi...	12324226	39	1411271
Caelo	Prednisolone, mi...	161343	44	20160601*
Caelo	Prednisolone, mi...	154000	60	20160225*
Caelo	Prednisolone, mi...	154000	59	20160225*
Caelo	Prednisolone, mi...	12324202	40	not necessary
Caelo	Prednisone	12232716	60	not necessary
Caelo	Prednisone	12232716	60	not necessary
Caelo	Prilocaine hydro...	15230305	60	not necessary
Caelo	Procaine hydroch...	15006906	60	not necessary
Caelo	Procaine hydroch...	13109906	60	AR-15-FG-008935-01
Caelo	Progesterone, mi...	17060903	40	not necessary
Caelo	Progesterone, mi...	15319601	40	not necessary
Caelo	Promethazine hyd...	161466	40	20160621*
Caelo	Promethazine hyd...	13119214	60	not necessary
Caelo	Propranolol hydr...	12370712	60	not necessary
Caelo	Propranolol hydr...	15283302	40	not necessary
Caelo	Propranolol hydr...	12370712	60	not necessary
Caelo	Propyl-4-hydroxy...	13365009	60	not necessary
Caelo	Propyl-4-hydroxy...	11373523	40	1402440
Caelo	Propyl-4-hydroxy...	162595	40	20161005*
Caelo	Propyl-4-hydroxy...	162595	40	20161005*
Caelo	Propyl-4-hydroxy...	190732	40	20190404*
Caelo	Propyphenazone	13332206	60	not necessary
Caelo	Propyphenazone	12039103	40	1402051
Caelo	Pyridoxine hydro...	162225	40	20160825*
Caelo	Pyridoxine hydro...	14262302	60	not necessary
Caelo	Quinine dihydroc...	172072	40	20170830*
Caelo	Quinine dihydroc...	172129	40	20173008*
Caelo	Quinine hydrochl...	16283302	40	1802260
Caelo	Quinine hydrochl...	15376201	60	not necessary
Caelo	Quinine sulfate ...	13024604	60	AR-16-FG-005472-01
Caelo	Quinine sulfate ...	15230005	40	not necessary
Caelo	Quinoline yellow	17034101	40	1709327
Caelo	Quinoline yellow	13109108	40	1711481
Caelo	Quinoline yellow	12042010	40	not necessary
Caelo	Quinoline yellow	172277	40	20171020*
Caelo	Quinolinol sulfa...	15147704	40	not necessary
Caelo	Quinolinol sulfa...	15147704	40	not necessary
Caelo	Quinolinol sulfa...	162218	40	20160808*
Caelo	Quinolinol sulfa...	13403501	60	not necessary
Caelo	Quinolinol sulfa...	13071104	40	1403095
Caelo	Resocrinol	13293108	60	AR-15-FG-006906-01
Caelo	Resocrinol	161533	40	20160524*
Caelo	Riboflavin	14322008	40	AR-16-FG-013669-01
Caelo	Riboflavin	13305203	60	AR-15-FG-005686-01
Caelo	Riboflavin	12003317	40	1405155

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Rice starch	16144501	40	not necessary
Caelo	Rice starch	12068525	40	not necessary
Caelo	Saccharin sodium	172592	40	20171206*
Caelo	Saccharin sodium	172110	40	20170831*
Caelo	Saccharin sodium	162210	39	20160811*
Caelo	Salicylic acid	13344018	60	not necessary
Caelo	Salicylic acid	13344002	60	not necessary
Caelo	Salicylic acid	160109	60	20160229*
Caelo	Selenium yeast	17002208	40	not necessary
Caelo	Selenium yeast	173332	40	20180207*
Caelo	Selenium yeast	17002202	40	1707512
Caelo	Selenium yeast	13347911	40	1509355
Caelo	Siam benzoin	12321312	40	pending
Caelo	Sildenafil citrate	170715	40	20170529*
Caelo	Sildenafil citrate	170715	40	20170529*
Caelo	Sildenafil citrate	160498	60	20160307*
Caelo	Sodium acetate t...	172003	40	20171107*
Caelo	Sodium acetate t...	12316506	40	1402111
Caelo	Sodium acetate t...	13411404	60	not necessary
Caelo	Sodium ascorbate	182620	40	20181112*
Caelo	Sodium ascorbate	15429905	40	not necessary
Caelo	Sodium benzoate	11267016	40	not necessary
Caelo	Sodium benzoate	13178619	60	not necessary
Caelo	Sodium citrate	13089335	60	AR-15-FG-000405-01
Caelo	Sodium citrate	12045026	40	not necessary
Caelo	Sodium citrate	14112517	60	not necessary
Caelo	Sodium cromoglic...	11114501	80	not necessary
Caelo	Sodium cromoglic...	15224404	60	not necessary
Caelo	Sodium cromoglic...	13232802	20	1505352
Caelo	Sodium cromoglic...	13232802	110	not necessary
Caelo	Sodium cyclamate	182167	40	20181023*
Caelo	Sodium cyclamate	12256803	40	1401418
Caelo	Sodium cyclamate	12105622	40	not necessary
Caelo	Sodium dodecyl s...	16229217	40	not necessary
Caelo	Sodium dodecyl s...	170955	40	20170503*
Caelo	Sodium dodecyl s...	15344503	60	not necessary
Caelo	Sodium dodecyl s...	190923	40	20190521*
Caelo	Sodium dodecyl s...	191243	40	20190521*
Caelo	Sodium dodecyl s...	13099011	40	1403108
Caelo	Sodium edetate	13142009	40	not necessary
Caelo	Sodium edetate	161364	60	20160504*
Caelo	Sodium methyl-4...	170519	40	20170308*
Caelo	Sodium methyl-4...	15358708	40	AR-16-FG-015095-01
Caelo	Sodium methyl-4...	15358713	40	not necessary
Caelo	Sodium methyl-4...	191439	40	20190628*
Caelo	Sodium methyl-4...	15358708	40	not necessary
Caelo	Sodium salicylat...	12158307	40	1403225
Caelo	Sodium salicylat...	13140513	60	not necessary
Caelo	Sodium salicylat...	163053	40	20161118*
Caelo	Sodium tetraborate	13290118	60	1412042
Caelo	Sodium tetraborate	161372	40	20160624*
Caelo	Sodium tetraborate	161372	60	20160624*
Caelo	Softisan® 601	15184905	60	not necessary
Caelo	Softisan® 601	13296104	40	1403144
Caelo	Softisan® 601	172178	40	20170817*
Caelo	Sorbic acid	15382002	40	not necessary
Caelo	Sorbic acid	172056	40	20170823*
Caelo	Sorbic acid	12038025	40	not necessary
Caelo	Sorbic acid	14176214	40	not necessary
Caelo	Sorbic acid	14380714	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Sorbic acid	191114	40	20190502*
Caelo	Sorbic acid	162080	40	20160722*
Caelo	Sorbic acid	13334008	60	AR-15-FG-007849-01
Caelo	Sorbic acid	162080	40	20160722*
Caelo	Sorbitan stearate	162352	40	20160826*
Caelo	Sorbitan stearate	17030805	40	not necessary
Caelo	Sorbitan stearate	172564	40	20171027*
Caelo	Sorbitan stearate	12274801	40	not necessary
Caelo	Sorbitol	13278217	40	not necessary
Caelo	Sorbitol	16038501	40	not necessary
Caelo	Sorbitol	16353107	40	not necessary
Caelo	Sorbitol	163531	40	20170119*
Caelo	Soy lecithin (in...	13279813	40	not necessary
Caelo	Soy lecithin (in...	17118101	40	not necessary
Caelo	Soy lecithin (in...	191879	40	20190712*
Caelo	Soy lecithin (in...	13030504	40	1404399
Caelo	Spironolactone	160773	60	20160415*
Caelo	Spironolactone	190158	40	20190308*
Caelo	Spironolactone	14378715	40	not necessary
Caelo	Spray-dried arab...	14141901	120	AR-15-FG-006809-01
Caelo	Spray-dried arab...	15413201	40	AR-17-FG-000844-01
Caelo	Spray-dried arab...	181860	40	20180725*
Caelo	Stearic acid 50	170810	40	20170426*
Caelo	Stearic acid 50	12241910	40	AR-15-FG-007284-01
Caelo	Stearic acid 50	14359314	60	not necessary
Caelo	Stearic acid 50	14359310	60	AR-16-FG-004044-01
Caelo	Stearic acid 50	14359310	60	not necessary
Caelo	Stearic acid 50	170810	40	20170426*
Caelo	Stemless carline...	12232802	40	not necessary
Caelo	Sucrose	173247	40	20180117*
Caelo	Sucrose	181582	40	20180605*
Caelo	Sucrose	13352909	40	1410501
Caelo	Sucrose	161999	40	20160714*
Caelo	Sucrose	171066	40	20170511*
Caelo	Sucrose	14340003	60	not necessary
Caelo	Sucrose	163286	40	20161219*
Caelo	Sucrose	162060	40	20160808*
Caelo	Sucrose	172133	40	20170817*
Caelo	Sucrose	163081	40	20161115*
Caelo	Sulfacetamide	15445901	40	not necessary
Caelo	Sulfacetamide	15445901	40	not necessary
Caelo	Sulfacetamide	11085910	40	1402651
Caelo	Sulfadiazine	17056602	40	not necessary
Caelo	Sulfadiazine	191850	40	20190730*
Caelo	Sulfadiazine	12164101	40	1402053
Caelo	Sulfadiazine	170566	40	20170330*
Caelo	Sulfadiazine	17056602	40	not necessary
Caelo	Sulfadiazine	13150011	60	not necessary
Caelo	Sulfogaiacol hyd...	190303	40	20190222*
Caelo	Sulfogaiacol hyd...	13405805	45	not necessary
Caelo	Tannin	12336016	60	not necessary
Caelo	Tannin	172092	40	20170829*
Caelo	Tartaric acid	161468	40	20160607*
Caelo	Tartaric acid	15130201	60	not necessary
Caelo	Tartaric acid	15164001	60	1704461
Caelo	Testosterone	13406004	50	not necessary
Caelo	Testosterone	15244906	38	not necessary
Caelo	Testosterone	16311506	40	not necessary
Caelo	Testosterone	13406004	50	not necessary
Caelo	Testosterone	15244906	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Testosterone pro...	12163607	40	AR-15-FG-000623-01
Caelo	Testosterone pro...	15444709	40	1707612
Caelo	Tetracaine base	181249	40	20180605*
Caelo	Tetracaine base	171117	40	20170512*
Caelo	Tetracaine hydro...	15437108	40	not necessary
Caelo	Tetracaine hydro...	12019822	40	not necessary
Caelo	Tetracycline hyd...	161167	60	20160504*
Caelo	Tetracycline hyd...	13201303	60	not necessary
Caelo	Tetracycline hyd...	13201303	60	not necessary
Caelo	Tetracycline hyd...	15230712	40	not necessary
Caelo	Thalidomide	13395704	60	not necessary
Caelo	Thalidomide	15276711	40	1704508
Caelo	Thalidomide	182496	40	20181022*
Caelo	Theophylline	12341101	60	not necessary
Caelo	Theophylline	12166316	40	1402055
Caelo	Theophylline-eth...	12145301	40	not necessary
Caelo	Theophylline-eth...	13268609	40	1412133
Caelo	Theophylline-eth...	12363708	40	not necessary
Caelo	Theophylline-eth...	17073707	30	not necessary
Caelo	Theophylline-eth...	15380302	40	not necessary
Caelo	Thiamine chlorid...	14069502	60	8986
Caelo	Thiamine chlorid...	14069502	60	not necessary
Caelo	Thiamine chlorid...	160156	60	20160316*
Caelo	Thymol	16144101	40	not necessary
Caelo	Thymol	14025805	60	not necessary
Caelo	Tosylchloramide ...	13103201	40	not necessary
Caelo	Tosylchloramide ...	13361602	40	1407634
Caelo	Tosylchloramide ...	16236603	40	not necessary
Caelo	Tretinoin	15239503	60	not necessary
Caelo	Tretinoin	15239503	60	not necessary
Caelo	Tretinoin	15028403	60	1507086
Caelo	Tretinoin	15028403	60	1507086
Caelo	Triamcinolone ac...	13379921	60	1503634
Caelo	Triamcinolone ac...	13256607	100	AR-15-FG-003755-01
Caelo	Triamcinolone, m...	15106905	45	not necessary
Caelo	Triamcinolone, m...	15106905	40	not necessary
Caelo	Triamcinolone, m...	15106911	40	not necessary
Caelo	Triclosan	14121209	45	not necessary
Caelo	Triclosan	13391808	60	not necessary
Caelo	Trometamol	14050002	40	not necessary
Caelo	Trometamol	16168810	40	not necessary
Caelo	Trometamol	161688	40	20160615*
Caelo	Ubidecarenone	172094	35	20170823*
Caelo	Ubidecarenone	14004504	60	AR-15-FG-009220-01
Caelo	Ubidecarenone	17209402	40	not necessary
Caelo	Ubidecarenone	17209402	40	not necessary
Caelo	Urea	18023003	40	not necessary
Caelo	Urea	161074	60	20160509*
Caelo	Urea	15154103	60	not necessary
Caelo	Urea	15154103	60	not necessary
Caelo	Urea	170216	40	20170202*
Caelo	Vanillin	160895	60	20160421*
Caelo	Vanillin	15135302	60	AR-17-FG-022168-03
Caelo	Vanillin	14044413	60	not necessary
Caelo	Vanillin	191201	40	20190515*
Caelo	Wax, bleached	13378619	20	1509364
Caelo	Wax, bleached	15292022	40	not necessary
Caelo	Wax, bleached	14336620	60	not necessary
Caelo	Wax, yellow	13335507	60	not necessary
Caelo	Wax, yellow	15317602	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Wheat starch	14023203	60	AR-16-FG-011433-01
Caelo	Wheat starch	15124707	60	not necessary
Caelo	Wool alcohols	191275	40	20190606*
Caelo	Wool alcohols	163530	40	20161223*
Caelo	Wool alcohols	13380601	40	1412037
Caelo	Wool alcohols	153703	60	20151111*
Caelo	Wool alcohols	153703	60	20151111*
Caelo	Xylitol	14068206	60	not necessary
Caelo	Xylitol	15135405	40	not necessary
Caelo	Xylitol	15135405	60	not necessary
Caelo	Xylometazoline h...	16074211	40	not necessary
Caelo	Xylometazoline h...	13344901	60	not necessary
Caelo	Xylometazoline h...	13344901	60	1501655
Caelo	Xylose	14218902	60	AR-14-FG-013086-01
Caelo	Xylose	17160902	40	not necessary
Caelo	Zinc acetate dihydroxide	170204	40	20170126*
Caelo	Zinc acetate dihydroxide	153833	40	20151127*
Caelo	Zinc acetate dihydroxide	15383301	40	not necessary
Caelo	Zinc gluconate	170951	40	20170518*
Caelo	Zinc gluconate	170165	40	20170125*
Caelo	Zinc stearate	180628	40	20180322*
Dr. Behr	L-carnitine hydrochloride	WE11030091	60	from supplier
Dr. Behr	L-carnitine tartaric acid	WE15080050	60	from supplier
Dr. Behr	L-carnitine tartaric acid	WE11060069	40	from supplier
Euro OTC	Acetylcysteine	1510016	60	20151116*
Euro OTC	Acetylcysteine	1208026-01	60	not necessary
Euro OTC	Acetylcysteine	1610011	40	20161107*
Euro OTC	Acetylsalicylic acid	1608037	40	20160921*
Euro OTC	Acetylsalicylic acid	1307024-01	45	not necessary
Euro OTC	Acetylsalicylic acid	1407026-01	60	not necessary
Euro OTC	Acetylsalicylic acid	1601042	40	20160204*
Euro OTC	Alfatradiol	1404001-01	50	1509288
Euro OTC	Alfatradiol	1603019	40	20160419*
Euro OTC	Alfatradiol	1603019	40	20160419*
Euro OTC	Alfatradiol	1404001-01	60	1509288
Euro OTC	Allantoin	1406011-01	60	not necessary
Euro OTC	Allantoin	1609012	40	20161006*
Euro OTC	Allantoin	1609012-01	40	not necessary
Euro OTC	Allantoin	1406011-04	40	not necessary
Euro OTC	alpha-tocopheryl acetate	1507001	60	20130708*
Euro OTC	alpha-tocopheryl acetate	1907033	40	20190820*
Euro OTC	alpha-tocopheryl acetate	1110021-02	40	not necessary
Euro OTC	alpha-tocopheryl acetate	1603046	40	20160429*
Euro OTC	alpha-tocopheryl acetate	1305032-01	60	AR-15-FG-001194-01
Euro OTC	Ambroxol hydrochloride	1604035	40	20160428*
Euro OTC	Ambroxol hydrochloride	1311032-01	40	not necessary
Euro OTC	Ambroxol hydrochloride	1602015	60	20160217*
Euro OTC	Ambroxol hydrochloride	11110008-01	40	1401146
Euro OTC	Ambroxol hydrochloride	1509032	60	20151027*
Euro OTC	Ambroxol hydrochloride	1601008	60	20160119*
Euro OTC	Amphotericin B	1512026-01	40	not necessary
Euro OTC	Amphotericin B	1305033-02	40	1403224
Euro OTC	Arginine	1702014	40	20170316*
Euro OTC	Arginine	1608031	40	20160916*
Euro OTC	Arginine	1604003	40	20160425*
Euro OTC	Arginine hydrochloride	1306023-01	60	AR-15-FG-008412-01
Euro OTC	Arginine hydrochloride	1904014	40	20190509*
Euro OTC	Arginine hydrochloride	1610042	40	20161128*
Euro OTC	Arginine hydrochloride	L1212009-01	40	1402430
Euro OTC	Arginine hydrochloride	L1103042-01	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Arginine hydroch...	1603033	40	20160425*
Euro OTC	Ascorbic acid	1607025	40	20160801*
Euro OTC	Ascorbic acid	1612018	40	20170109*
Euro OTC	Ascorbic acid	1605024	40	20160617*
Euro OTC	Ascorbic acid	1609015	40	20161006*
Euro OTC	Beclometasone di...	1203022-03	60	not necessary
Euro OTC	Beclometasone di...	1502032	40	20150327*
Euro OTC	Beclometasone di...	1502032	30	20150327*
Euro OTC	Beclometasone di...	1203022-03	55	1509281
Euro OTC	Betamethasone di...	1601034	45	20160204*
Euro OTC	Betamethasone di...	1311023-01	60	1505455
Euro OTC	Betamethasone di...	1601034	45	20160204*
Euro OTC	Betamethasone di...	1311023-01	46	1505455
Euro OTC	Betamethasone va...	1811008	30	20181207*
Euro OTC	Betamethasone va...	1510037	60	20151215*
Euro OTC	Betamethasone va...	1610010	40	20161116*
Euro OTC	Betamethasone va...	1406019-01	60	not necessary
Euro OTC	Betamethasone va...	1811008	40	20181207*
Euro OTC	Betamethasone va...	1208014-01	40	not necessary
Euro OTC	Betamethasone va...	1510037	30	20151215*
Euro OTC	Bifonazole	1401031-01	60	not necessary
Euro OTC	Bifonazole	1510001	60	20151203*
Euro OTC	Bifonazole	1206027-01	40	1312397
Euro OTC	Biotin	1402028-01	60	not necessary
Euro OTC	Biotin	1610003	40	20161110*
Euro OTC	Biotin	1508030	60	20150925*
Euro OTC	Budesonide, micr...	1311022-01	60	not necessary
Euro OTC	Budesonide, micr...	1604018	23	20160509*
Euro OTC	Budesonide, micr...	1311022-01	37	not necessary
Euro OTC	Budesonide, micr...	1604018	40	20160509*
Euro OTC	Caffeine	1808017	40	20180905*
Euro OTC	Caffeine	1211012-01	40	1412248
Euro OTC	Calcium ascorbate	1709004	40	20170914*
Euro OTC	Calcium ascorbate	1709005	40	20170914*
Euro OTC	Calcium ascorbate	1801010	40	20180115*
Euro OTC	Calcium ascorbate	1612003	40	20170110*
Euro OTC	Calcium citrate ...	1509039	60	20151013*
Euro OTC	Calcium citrate ...	1701001	40	20170124*
Euro OTC	Calcium citrate ...	1509039-01	60	not necessary
Euro OTC	Calcium citrate ...	1305018-01	40	not necessary
Euro OTC	Calcium gluconate	1407041-01	60	not necessary
Euro OTC	Calcium gluconate	1304012-01	40	AR-15-FG-008127-01
Euro OTC	Calcium orotate ...	1603003	40	20160329*
Euro OTC	Calcium orotate ...	1401006-01	40	1412540
Euro OTC	Calcium orotate ...	1603003	40	20180615*
Euro OTC	Calcium pantothe...	1305034-01	40	not necessary
Euro OTC	Calcium pantothe...	1607023	40	20160809*
Euro OTC	Calcium pantothe...	1509024	60	20151005*
Euro OTC	Calcium pantothe...	1702018	40	20170330*
Euro OTC	Calcium pantothe...	1711003	40	20171123*
Euro OTC	Calcium pantothe...	1603042	40	20160429*
Euro OTC	Camphor, racemic	1406029-01	60	not necessary
Euro OTC	Capsaicin	1811009	40	20181207*
Euro OTC	Capsaicin	1703023	40	20170407*
Euro OTC	Carbomers 35000 ...	1610024	40	20161102*
Euro OTC	Carbomers 35000 ...	1607028	40	20160804*
Euro OTC	Carbomers 35000 ...	1604032	40	20160428*
Euro OTC	Carbomers 35000 ...	1512010	60	20151216*
Euro OTC	Carbomers 35000 ...	1311027-02	60	not necessary
Euro OTC	Cellulose, micro...	1211037-01	30	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Cellulose, micro...	1211037-01	40	not necessary
Euro OTC	Cellulose, micro...	1201004-01	40	not necessary
Euro OTC	Cellulose, micro...	1507022	60	20150805*
Euro OTC	Cellulose, micro...	1604029	40	20160512*
Euro OTC	Cetyl alcohol	1508014-01	40	not necessary
Euro OTC	Chloramphenicol	1407043-01	60	not necessary
Euro OTC	Chloramphenicol	1209019-02	80	not necessary
Euro OTC	Chloramphenicol	1209019-02	80	not necessary
Euro OTC	Chloramphenicol	1602002	45	20160309*
Euro OTC	Chlortetracyclin...	1604015	40	20160510*
Euro OTC	Chlortetracyclin...	1108051-02	40	not necessary
Euro OTC	Cholecalciferol ...	1608002	40	20160824*
Euro OTC	Cholecalciferol ...	1602019	40	20160307*
Euro OTC	Cholecalciferol ...	1701053	40	20170208*
Euro OTC	Cholecalciferol ...	1801038	40	20180206*
Euro OTC	Cholecalciferol ...	1802014	40	20180223*
Euro OTC	Cholecalciferol ...	1502021-01	60	not necessary
Euro OTC	Cholecalciferol ...	1707027	40	20171020*
Euro OTC	Cholecalciferol ...	1304016-02	40	not necessary
Euro OTC	Cholecalciferol ...	1906009	40	20190619*
Euro OTC	Cholecalciferol ...	1210051-01	40	1403341
Euro OTC	Cholecalciferol ...	1508032	60	20150901*
Euro OTC	Choline hydrogen...	1603044-01	40	not necessary
Euro OTC	Choline hydrogen...	1901053	40	20190228*
Euro OTC	Choline hydrogen...	1702005	40	20170307*
Euro OTC	Choline hydrogen...	1607020	40	20160808*
Euro OTC	Choline hydrogen...	1603044-01	40	not necessary
Euro OTC	Choline hydrogen...	1603044	40	20160429*
Euro OTC	Citrulline	1604038	40	20160523*
Euro OTC	Citrulline	1703036	40	20170418*
Euro OTC	Citrulline	1403021-01	59	1411031
Euro OTC	Citrulline	1811028	40	20181127*
Euro OTC	Citrulline	1310015-01	40	not necessary
Euro OTC	Citrulline	1506014	60	20150624*
Euro OTC	Clindamycin hydr...	1511004	60	20151123*
Euro OTC	Clindamycin hydr...	1510038	60	20151120*
Euro OTC	Clindamycin hydr...	1111058-01	80	not necessary
Euro OTC	Clioquinol	1309009-02	40	AR-15-FG-009542-01
Euro OTC	Clioquinol	1509038	60	20151009*
Euro OTC	Clobetasol propi...	1211021-04	60	1605044
Euro OTC	Clobetasol propi...	1606022	40	20160713*
Euro OTC	Clobetasol propi...	1606022	40	20160713*
Euro OTC	Clotrimazole	1510003	60	20151111*
Euro OTC	Clotrimazole	1306036-01	40	not necessary
Euro OTC	Colistin sulfate	1507044	60	20150817*
Euro OTC	Colistin sulfate	1602021	60	20160307*
Euro OTC	Colistin sulfate	1406001-01	60	1508600
Euro OTC	Colistin sulfate	1606008	40	20160704*
Euro OTC	Cyanocobalamin	1812011	40	20190205*
Euro OTC	Cyanocobalamin	1604017	30	20160509*
Euro OTC	Cyanocobalamin	1604017	30	20160509*
Euro OTC	Cyanocobalamin	1506035	60	0000010479
Euro OTC	Cyanocobalamin	1506035	60	20150708*
Euro OTC	Cyanocobalamin	1506035	45	0000010479
Euro OTC	Cyanocobalamin	1506035	60	20150708*
Euro OTC	Cyanocobalamin	1305040-01	60	AR-15-FG-005214-01
Euro OTC	Cyproterone acet...	1611003	40	20161201*
Euro OTC	Cyproterone acet...	1810017	40	20181108*
Euro OTC	Cystine	1610013	40	20161110*
Euro OTC	Cystine	1408016-01	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Cystine	1307052-01	40	not necessary
Euro OTC	Cystine	1603010	40	20160406*
Euro OTC	Cystine	1801016	40	20180126*
Euro OTC	Deoxycholic acid	1305051-01	40	not necessary
Euro OTC	Dequalinium chlo...	1405012-01	60	not necessary
Euro OTC	Dequalinium chlo...	1512032	60	20160121*
Euro OTC	Dequalinium chlo...	1512032-02	40	1809202
Euro OTC	Dexamethasone	1611001	40	20161212*
Euro OTC	Dexamethasone	1601052	45	20160225*
Euro OTC	Dexamethasone	1303007-01	60	not necessary
Euro OTC	Dexamethasone	1601052	45	20160225*
Euro OTC	Dexamethasone	1303007-01	40	1407432
Euro OTC	Dexamethasone	1611001	37	20161212*
Euro OTC	Dexamethasone ac...	1405024-01	45	not necessary
Euro OTC	Dexamethasone ac...	1708011	40	20171002*
Euro OTC	Dexamethasone ac...	1405024-01	60	not necessary
Euro OTC	Diclofenac sodium	1607016	40	20160808*
Euro OTC	Diclofenac sodium	1403013-01	40	not necessary
Euro OTC	Dimethyl fumarate	1410013-01	60	not necessary
Euro OTC	Dithranol	1702022-01	40	1809448
Euro OTC	EDTA calcium dis...	1705029	40	20170626*
Euro OTC	EDTA calcium dis...	1807024	40	20180725*
Euro OTC	Emulsifying ceto...	1510043	60	20151124*
Euro OTC	Emulsifying ceto...	1504010-01	60	not necessary
Euro OTC	Erythromycin	1308016-02	60	1412542
Euro OTC	Erythromycin	1705038	40	20170705*
Euro OTC	Erythromycin	1508034	60	20150925*
Euro OTC	Erythromycin	1607038	40	20160831*
Euro OTC	Erythromycin	1811019	40	20190114*
Euro OTC	Escin soluble	1605021	40	20160609*
Euro OTC	Escin soluble	1310020-01	40	AR-15-FG-005777-01
Euro OTC	Escin soluble	1506038	60	20150708*
Euro OTC	Estradiol benzoate	1507009	60	0000010566
Euro OTC	Estradiol benzoate	1507009	60	0000010566
Euro OTC	Estradiol benzoate	1301006-02	40	1507405
Euro OTC	Estradiol hemihy...	1507013	60	1510623
Euro OTC	Estradiol hemihy...	1507013	60	1510623
Euro OTC	Estradiol hemihy...	1612008	40	20170117*
Euro OTC	Estradiol hemihy...	1510039	89	20151103*
Euro OTC	Estradiol hemihy...	1208022-01	40	not necessary
Euro OTC	Estradiol hemihy...	1612008	40	20170117*
Euro OTC	Estradiol hemihy...	1510039	60	20151103*
Euro OTC	Estradiol hemihy...	1601017	45	20160122*
Euro OTC	Estriol	1509042	60	20151009*
Euro OTC	Estriol	1507026	60	0000010671
Euro OTC	Estriol	1507026	60	0000010671
Euro OTC	Estriol	1703044	40	20170502*
Euro OTC	Estriol	1509042	60	20151009*
Euro OTC	Ethacridine lact...	1412034-02	40	not necessary
Euro OTC	Ethacridine lact...	1206016-02	60	not necessary
Euro OTC	Ethyl hydrogen f...	1510005	60	20151030*
Euro OTC	Ethyl hydrogen f...	1310029-01	60	not necessary
Euro OTC	Ethyl hydrogen f...	1306012-01	60	1505456
Euro OTC	Ethyl hydrogen f...	1510006-01	40	AR-18-FG-000547-01
Euro OTC	Ethyl hydrogen f...	1601010	60	20160114*
Euro OTC	Ethyl hydrogen f...	1510006	60	20151030*
Euro OTC	Finasteride (for...	1605020	40	20160613*
Euro OTC	Finasteride (for...	1605020	40	20160613*
Euro OTC	Finasteride (for...	1403016-01	60	1506004
Euro OTC	Finasteride (for...	1403016-02	51	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Finasteride (for...	1610012	33	20161116*
Euro OTC	Finasteride (for...	1610012	40	20161116*
Euro OTC	Finasteride (for...	1403016-02	60	not necessary
Euro OTC	Fluocinolonaceto...	2804013E-01	60	not necessary
Euro OTC	Fluocinolonaceto...	2804013E-01	60	1404008
Euro OTC	Fluocinolonaceto...	2804013E-01	60	not necessary
Euro OTC	Fluocinolonaceto...	2804013E-01	60	not necessary
Euro OTC	Folic acid hydrate	1508024	60	20150902*
Euro OTC	Folic acid hydrate	1604036	40	20160429*
Euro OTC	Folic acid hydrate	1611009	40	20161201*
Euro OTC	Folic acid hydrate	1401005-01	40	not necessary
Euro OTC	Fructose	1509037	60	20151027*
Euro OTC	Fructose	1212015-01	40	not necessary
Euro OTC	Fructose	1512009	60	20160113*
Euro OTC	Fructose	1701049	40	20170213*
Euro OTC	Gentamicin sulfate	1506016	60	0000010321
Euro OTC	Gentamicin sulfate	1711025	40	20171206*
Euro OTC	Gentamicin sulfate	1305047-03	60	1411076
Euro OTC	Gentamicin sulfate	1305047-03	45	AR-15-FG-012403-01
Euro OTC	Gentamicin sulfate	1907036	40	20190823*
Euro OTC	Gentamicin sulfate	1305047-01	30	not necessary
Euro OTC	Gentamicin sulfate	1610036	40	20161124*
Euro OTC	Gentamicin sulfate	1610036-01	40	1804343
Euro OTC	Glucose monohydr...	1807005	40	20180718*
Euro OTC	Glucose monohydr...	1701002	40	20170213*
Euro OTC	Glucose monohydr...	1509006	60	20150914*
Euro OTC	Glucose monohydr...	1604006	40	20160511*
Euro OTC	Glucose monohydr...	1601016	60	20160129*
Euro OTC	Glutamic acid	1604056	40	20180613*
Euro OTC	Glutamic acid	1604056	40	20160523*
Euro OTC	Glutamic acid	L1405023-01	60	not necessary
Euro OTC	Glutamine	1706010	40	20170802*
Euro OTC	Glutamine	1702007	40	20170322*
Euro OTC	Glycine	1801005	40	20180125*
Euro OTC	Glycine	1402022-01	60	AR-15-FG-005430-01
Euro OTC	Glycine	1402022-01	40	not necessary
Euro OTC	Griffonia seed e...	1905017	40	20190709*
Euro OTC	Griffonia seed e...	1205027	40	from supplier
Euro OTC	Griffonia seed e...	1705007	40	20170601*
Euro OTC	Guaifenesin	1212002-01	40	not necessary
Euro OTC	Guaifenesin	1506044	60	20150708*
Euro OTC	Guaifenesin	1605005	40	20160525*
Euro OTC	Heparin sodium	12050009-03	40	AR-15-FG-007427-02
Euro OTC	Heparin sodium	1605022-01	40	9699
Euro OTC	Histamine dihydr...	1111004-02	40	1403674
Euro OTC	Histamine dihydr...	1506030	60	20150706*
Euro OTC	Histidine	1703022	40	20170405*
Euro OTC	Hydrocortisone	1206018-01	40	not necessary
Euro OTC	Hydrocortisone	1509029	60	20151014*
Euro OTC	Hydrocortisone	1509029	51	20151014*
Euro OTC	Hydrocortisone	1611016	40	20161220*
Euro OTC	Hydrocortisone a...	1508029	60	20150902*
Euro OTC	Hydrocortisone a...	1604052	40	20160525*
Euro OTC	Hydrocortisone a...	1209003-01	40	not necessary
Euro OTC	Hydrocortisone a...	1209003-01	45	not necessary
Euro OTC	Hydrocortisone a...	1508029	44	20150902*
Euro OTC	Hydroxyethylcell...	1110024-04	40	not necessary
Euro OTC	Hydroxyethylcell...	1804003	40	20180413*
Euro OTC	Ibuprofen	1211011-02	40	1501533
Euro OTC	Ibuprofen	1410025-02	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Indometacin	1205002-01	40	not necessary
Euro OTC	Indometacin	1512031	60	20160122*
Euro OTC	Indometacin	1612017	40	20170209*
Euro OTC	Indometacin	1012009-02	40	1402050
Euro OTC	Inositol	1206032-01	40	not necessary
Euro OTC	Isoxuprine hydro...	1503008-01	40	not necessary
Euro OTC	Lactose monohydr...	1509030	60	20151013*
Euro OTC	Lactose monohydr...	1606025	40	20160628*
Euro OTC	L-carnitine tart...	1302017-01	40	from supplier
Euro OTC	Levocarnitine	1210029-01	40	not necessary
Euro OTC	L-glutathione	1808036	40	20180919*
Euro OTC	L-glutathione	1703008	40	20170327*
Euro OTC	L-glutathione	1102019-02	60	AR-15-FG-008218-01
Euro OTC	Lidocaine	1612021	40	20170206*
Euro OTC	Lidocaine	1310023-01	60	not necessary
Euro OTC	Lidocaine hydroc...	1907012	40	20190801*
Euro OTC	Lidocaine hydroc...	1403025-02	60	not necessary
Euro OTC	Lidocaine hydroc...	1411010	60	20150126*
Euro OTC	L-isoleucine	1704022	40	20170612*
Euro OTC	L-isoleucine	1604057	40	20160601*
Euro OTC	L-leucine	1302013-01	60	not necessary
Euro OTC	L-leucine	1601032	40	20160308*
Euro OTC	L-phenylalanine	1704003	40	20170425*
Euro OTC	L-phenylalanine	1212014-03	60	not necessary
Euro OTC	L-proline	1606010	40	20160714*
Euro OTC	L-proline	1705021	40	20170720*
Euro OTC	L-threonine	1210027-01	60	not necessary
Euro OTC	L-threonine	1210027-01	40	1403152
Euro OTC	L-threonine	1601031	40	20160225*
Euro OTC	L-threonine	1909003	40	20190913*
Euro OTC	L-tryptophan	1506003	60	20150630*
Euro OTC	L-tryptophan	1712006	40	20180105*
Euro OTC	L-tryptophan	1506003	60	20150630*
Euro OTC	Lysine hydrochlo...	1611021	40	20161202*
Euro OTC	Lysine hydrochlo...	1708009	40	20170922*
Euro OTC	Lysine hydrochlo...	1805030	40	20180611*
Euro OTC	Macrogols 1500	1508002	60	20150810*
Euro OTC	Macrogols 1500	1604005	40	20160512*
Euro OTC	Macrogols 4000	1605025	40	20160603*
Euro OTC	Macrogols 4000	1507004	60	0000010550
Euro OTC	Macrogols 4000	1711016	40	20171214*
Euro OTC	Macrogols 4000	1508031	60	20150914*
Euro OTC	Macrogols 4000	1702003	40	20170301*
Euro OTC	Macrogols 4000	1605012	40	20160606*
Euro OTC	Magnesium citrate	1602028	40	20160314*
Euro OTC	Magnesium citrate	1808020	40	20180827*
Euro OTC	Magnesium citrate	1511018	60	20151130*
Euro OTC	Magnesium citrate	1711018	40	20171129*
Euro OTC	Magnesium citrate	1506020	60	20150623*
Euro OTC	Magnesium citrate	1501031	60	AR-15-FG-007845-01
Euro OTC	Magnesium citrate	1706007	40	20170720*
Euro OTC	Magnesium citrate	1902006	40	20190226*
Euro OTC	Magnesium citrate	1803033	40	20180403*
Euro OTC	Magnesium citrate	1507019	60	20150722*
Euro OTC	Magnesium hydrog...	1906029	40	20190712*
Euro OTC	Magnesium hydrog...	1901036	40	20190227*
Euro OTC	Magnesium hydrog...	1204014-01	40	1404588
Euro OTC	Magnesium hydrog...	1306045-01	60	not necessary
Euro OTC	Magnesium hydrog...	1605031	40	20160617*
Euro OTC	Magnesium hydrog...	1601005	60	20160119*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Magnesium hydrog...	1510033	60	20151117*
Euro OTC	Magnesium hydrog...	1507010	60	20150713*
Euro OTC	Magnesium hydrog...	1701003	40	20170124*
Euro OTC	Magnesium orotat...	1511021	60	20151228*
Euro OTC	Magnesium orotat...	1707018	40	20170907*
Euro OTC	Magnesium stearate	1502033-01	40	not necessary
Euro OTC	Magnesium stearate	1704002	40	20170622*
Euro OTC	Manganese orotat...	1907009	40	20190718*
Euro OTC	Manganese orotat...	1601003	60	20160113*
Euro OTC	Manganese orotat...	1801031	40	20180131*
Euro OTC	Mannitol	1508012	60	20150827*
Euro OTC	Mannitol	1511025	60	20151228*
Euro OTC	Melatonin	1709020	40	20180105*
Euro OTC	Melatonin	1804015	40	20180522*
Euro OTC	Melatonin	1804015	40	20180522*
Euro OTC	Melatonin	1510041	60	20151123*
Euro OTC	Melatonin	1305029-02	60	1411054
Euro OTC	Menthol	1403022-02	45	not necessary
Euro OTC	Menthol	1608012	40	20160902*
Euro OTC	Menthol	1603024	40	20160429*
Euro OTC	Menthol	1509045	60	20151008*
Euro OTC	Mesalazine	1206029-04	60	not necessary
Euro OTC	Mesalazine	1206029-04	60	AR-15-FG-011819-01
Euro OTC	Mesalazine	1602005-01	40	not necessary
Euro OTC	Mesalazine	1602005	40	20160218*
Euro OTC	Methionine	1602012	60	20160219*
Euro OTC	Methionine	L1404041-01	60	not necessary
Euro OTC	Methoxsalen	1506013	60	20150706*
Euro OTC	Methoxsalen	130604101	60	not necessary
Euro OTC	Methoxsalen	1608019	40	20160913*
Euro OTC	Methyl-4-hydroxy...	1601037	60	20160119*
Euro OTC	Methyl-4-hydroxy...	1205025-01	40	not necessary
Euro OTC	Methylcobalamin	1601043	60	20160209*
Euro OTC	Methylcobalamin	1510002	60	20151201*
Euro OTC	Methylcobalamin	1207019-01	40	1402501
Euro OTC	Methylcobalamin	1809001	40	20180925*
Euro OTC	Methylcobalamin	1401032-01	60	8778
Euro OTC	Methylcobalamin	1804006	40	20180420*
Euro OTC	Methylcobalamin	1401032	60	20151201*
Euro OTC	Methylcobalamin	1904002	40	20190806*
Euro OTC	Metronidazole	1212030-02	40	not necessary
Euro OTC	Metronidazole	1508004	60	20150819*
Euro OTC	Miconazole nitrate	1307025-02	60	not necessary
Euro OTC	Miconazole nitrate	1509014	60	20150930*
Euro OTC	Minoxidil	1509035	60	20151016*
Euro OTC	Minoxidil	1612011	40	20170130*
Euro OTC	Minoxidil	1602013	60	20160225*
Euro OTC	Minoxidil	1311038-01	59	AR-15-FG-007848-01
Euro OTC	Minoxidil	1411006-01	60	not necessary
Euro OTC	Mometasone furoa...	1510036	60	20151118*
Euro OTC	Mometasone furoa...	1403015-01	31	not necessary
Euro OTC	Mometasone furoa...	1403015-01	60	not necessary
Euro OTC	Mometasone furoa...	1510036	30	20151118*
Euro OTC	Mometasone furoa...	1610009	40	20170321*
Euro OTC	Naphazoline hydr...	1506015	60	20150727*
Euro OTC	Naphazoline hydr...	1407034-01	60	AR-15-FG-007847-01
Euro OTC	Naphazoline nitr...	1810016	40	20181112*
Euro OTC	Naphazoline nitr...	1512033	60	20160126*
Euro OTC	Neomycin sulfate	1303033-03	60	AR-15-FG-016582-01
Euro OTC	Neomycin sulfate	1601046	45	20160212*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Nicotinamide	1511010	60	20151203*
Euro OTC	Nicotinamide	1601045	60	20160210*
Euro OTC	Nicotinamide	1610019	40	20161202*
Euro OTC	Nicotinamide	1310025-01	60	not necessary
Euro OTC	Nicotinamide	1906027	40	20190712*
Euro OTC	Nicotinic acid	1906026	40	20190711*
Euro OTC	Norethisterone a...	1306033-01	60	1511684
Euro OTC	Nystatin	1905014	40	20190709*
Euro OTC	Nystatin	1712010-01	40	1807445
Euro OTC	Ornithine aspart...	L1304034-01	40	1402571
Euro OTC	Ornithine aspart...	1611018-01	40	not necessary
Euro OTC	Ornithine aspart...	1611018	40	20161213*
Euro OTC	Ornithine aspart...	1312030-01	60	1506280
Euro OTC	Ornithine hydroc...	1307028-01	120	not necessary
Euro OTC	Ornithine hydroc...	1605011	40	20160601*
Euro OTC	Oxytetracycline ...	1506037	60	20150721*
Euro OTC	Oxytetracycline ...	1602025	40	20160609*
Euro OTC	Oxytetracycline ...	1701034	40	20170314*
Euro OTC	Oxytetracycline ...	1701035-02	40	not necessary
Euro OTC	Oxytetracycline ...	1306022-02	60	1507542
Euro OTC	Oxytetracycline ...	1506037	60	20150721*
Euro OTC	Oxytetracycline ...	1602025	40	20160609*
Euro OTC	Oxytetracycline ...	1602025-02	40	not necessary
Euro OTC	Oxytetracycline ...	1602025-01	40	AR-16-FG-014422-01
Euro OTC	Oxytetracycline ...	1602025-02	40	not necessary
Euro OTC	Paracetamol	1111024-02	60	not necessary
Euro OTC	Phenylbutazon	1507023	60	20150805*
Euro OTC	Phenylbutazon	1305050-02	60	AR-15-FG-005687-01
Euro OTC	Phenylephrine hy...	1704009	40	20170524*
Euro OTC	Pilocarpine hydr...	1305044-01	60	not necessary
Euro OTC	Pilocarpine hydr...	1603047	40	20160429*
Euro OTC	Pilocarpine hydr...	1305044-01	60	not necessary
Euro OTC	Pilocarpine hydr...	1603047	40	20160429*
Euro OTC	Polymyxin B sulf...	1304033-01	40	not necessary
Euro OTC	Polymyxin B sulf...	1505005-01	60	AR-16-FG-005387-01
Euro OTC	Polymyxin B sulf...	1604016	40	20160509*
Euro OTC	Polymyxin B sulf...	1202039-02	40	1411033
Euro OTC	Polymyxin B sulf...	1304033-02	60	not necessary
Euro OTC	Potassium citrate	1509044	60	20151013*
Euro OTC	Potassium citrate	1611017	40	20161207*
Euro OTC	Potassium sorbate	1907023	40	20190802*
Euro OTC	Potassium sorbate	1608013	40	20160819*
Euro OTC	Potassium sorbate	1508025	60	20150902*
Euro OTC	Povidone iodine	1606026-02	40	not necessary
Euro OTC	Povidone iodine	1509009	60	20150918*
Euro OTC	Povidone iodine	1501029-01	60	not necessary
Euro OTC	Povidone iodine	1905008	40	20190628*
Euro OTC	Povidone iodine	1805007	40	20180525*
Euro OTC	Prasterone	1102036-01	40	1404310
Euro OTC	Prasterone	1712014	40	20180115*
Euro OTC	Prasterone	1704012	40	20170529*
Euro OTC	Prasterone	1710035	40	20171212*
Euro OTC	Prasterone	1602018	30	20160301*
Euro OTC	Prasterone	1506011	60	20150629*
Euro OTC	Prasterone	1707014	40	20170829*
Euro OTC	Praziquantel	1908004	40	20190830*
Euro OTC	Prednicarbate, m...	1302033-02	60	1509286
Euro OTC	Prednicarbate, m...	1601058	45	20160223*
Euro OTC	Prednicarbate, m...	1507011	60	20150730*
Euro OTC	Prednicarbate, m...	1601058	31	20160223*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Prednisolone ace...	1512020	45	20160119*
Euro OTC	Prednisolone ace...	1512020	44	20160119*
Euro OTC	Prednisolone ace...	1306040	40	not necessary
Euro OTC	Prednisolone ace...	1209028-03	60	AR-15-FG-010480-01
Euro OTC	Prednisolone, mi...	1607011	40	20160818*
Euro OTC	Prednisolone, mi...	1212021-02	60	1411270
Euro OTC	Prednisolone, mi...	1212021-02	48	1411270
Euro OTC	Prednisolone, mi...	1607011	40	20160818*
Euro OTC	Prednisolone, mi...	1501047-01	30	AR-15-FG-015509-01
Euro OTC	Prednisolone, mi...	1501047-01	60	AR-15-FG-015509-01
Euro OTC	Prednisone	1401042-02	60	1605276
Euro OTC	Prednisone	1604051	40	20160523*
Euro OTC	Prednisone	1604051	40	20160523*
Euro OTC	Pregnenolone	1612019	40	20170131*
Euro OTC	Pregnenolone	1602001	45	20160208*
Euro OTC	Pregnenolone	1204040-01	40	1404311
Euro OTC	Pregnenolone	1312007-01	60	8985
Euro OTC	Procaine hydroch...	1601014	60	20160203*
Euro OTC	Progesterone, mi...	1506033	60	0000010482
Euro OTC	Progesterone, mi...	1506034	60	0000010483
Euro OTC	Progesterone, mi...	1511008	60	20151119*
Euro OTC	Progesterone, mi...	1601001	45	20160128*
Euro OTC	Progesterone, mi...	1601009	45	20160128*
Euro OTC	Progesterone, mi...	1601009	45	20160128*
Euro OTC	Progesterone, mi...	1603004	40	20160406*
Euro OTC	Progesterone, mi...	1509008	60	20150914*
Euro OTC	Progesterone, mi...	1506033	30	0000010482
Euro OTC	Progesterone, mi...	1111053-01	30	not necessary
Euro OTC	Pyridoxal-5-phos...	1701015-01	40	not necessary
Euro OTC	Pyridoxal-5-phos...	1902022	40	20190320*
Euro OTC	Pyridoxal-5-phos...	1701015	40	20170220*
Euro OTC	Pyridoxal-5-phos...	1302029-01	40	not necessary
Euro OTC	Pyridoxal-5-phos...	1908001	40	20190903*
Euro OTC	Pyridoxal-5-phos...	1302029-01	40	8897
Euro OTC	Pyridoxine hydro...	1407009-01	60	not necessary
Euro OTC	Pyridoxine hydro...	1601039	60	20160203*
Euro OTC	Pyridoxine hydro...	1507012	60	20150721*
Euro OTC	Pyridoxine hydro...	1401044-01	60	not necessary
Euro OTC	Ribavirin	1604012	40	20160428*
Euro OTC	Ribavirin	1308001	45	20130902*
Euro OTC	Riboflavin	1701050	40	20170214*
Euro OTC	Riboflavin	1701050-01	40	not necessary
Euro OTC	Salicylic acid	1110020-01	40	not necessary
Euro OTC	Salicylic acid	1509007	60	20150929*
Euro OTC	Selenium yeast	L1204012-01	40	1403461
Euro OTC	Serine	1512012	60	20160112*
Euro OTC	Serine	1606020	40	20160714*
Euro OTC	Serine	1606020-01	40	not necessary
Euro OTC	Serine	1711033	40	20171208*
Euro OTC	Serine	1904034	40	20190521*
Euro OTC	Sildenafil citrate	1311010-04	40	not necessary
Euro OTC	Sildenafil citrate	1311010-01	60	not necessary
Euro OTC	Sildenafil citrate	1908006	40	20190906*
Euro OTC	Sodium acetate t...	1706008	40	20170710*
Euro OTC	Sodium acetate t...	1507020	60	20150723*
Euro OTC	Sodium ascorbate	1406018-01	60	1411165
Euro OTC	Sodium ascorbate	1506018	60	20150706*
Euro OTC	Sodium ascorbate	1507030	60	20150811*
Euro OTC	Sodium ascorbate	1510017	60	20151102*
Euro OTC	Sodium ascorbate	1406018-01	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Sodium ascorbate	1601007	60	20160120*
Euro OTC	Sodium benzoate	1510044	60	20151109*
Euro OTC	Sodium benzoate	1409040-01	60	not necessary
Euro OTC	Sodium benzoate	1610037	40	20161202*
Euro OTC	Sodium benzoate	1506019	60	20150623*
Euro OTC	Sodium benzoate	1602020	60	20160225*
Euro OTC	Sodium benzoate	1701045	40	20170207*
Euro OTC	Sodium benzoate	1711019	40	20171205*
Euro OTC	Sodium citrate	1701029	40	20170126*
Euro OTC	Sodium citrate	1610016	40	20161102*
Euro OTC	Sodium citrate	1506041	60	0000010492
Euro OTC	Sodium citrate	1511017	60	20151130*
Euro OTC	Sodium cromoglic...	1703027	40	20170405*
Euro OTC	Sodium cromoglic...	1403017-01	60	not necessary
Euro OTC	Sodium cromoglic...	1510018	60	20151109*
Euro OTC	Sodium edetate	1812016	40	20190116*
Euro OTC	Sodium edetate	1609023	40	20170106*
Euro OTC	Sodium tetraborate	1404025-01	60	not necessary
Euro OTC	Sodium tetraborate	1510032	60	20151118*
Euro OTC	Sodium tetraborate	1711002	40	20171128*
Euro OTC	Sodium tetraborate	1507024	60	20150730*
Euro OTC	Sorbitol	1610018	40	20161121*
Euro OTC	Sorbitol	1506052	60	20150713*
Euro OTC	Sorbitol	1604007	40	20160428*
Euro OTC	Sorbitol	1512016	60	20160113*
Euro OTC	Spironolactone	1212034-01	42	not necessary
Euro OTC	Spironolactone	1602006	45	20160309*
Euro OTC	Spironolactone	1212034-01	50	1403151
Euro OTC	Sulfacetamide	1312012-01	60	not necessary
Euro OTC	Sulfacetamide	1509018	60	20150929*
Euro OTC	Tartaric acid	1907022	40	20190813*
Euro OTC	Taurine	1704001	40	20170426*
Euro OTC	Taurine	1607035	40	20160819*
Euro OTC	Taurine	1306024-01	50	not necessary
Euro OTC	Taurine	1904003	40	20190425*
Euro OTC	Taurine	1504032-01	60	9167
Euro OTC	Taurine	1602011	60	20160218*
Euro OTC	Taurine	1504032-01	60	not necessary
Euro OTC	Testosterone	1604049	34	20160510*
Euro OTC	Testosterone	1511005	35	20151119*
Euro OTC	Testosterone	1307001-01	40	1402309
Euro OTC	Testosterone	1604049	40	20160510*
Euro OTC	Testosterone	1511005	60	20151119*
Euro OTC	Testosterone pro...	1511003	43	20151111*
Euro OTC	Testosterone pro...	1511003	60	20151111*
Euro OTC	Testosterone pro...	1512023	45	20151230*
Euro OTC	Testosterone pro...	1512023	45	20151230*
Euro OTC	Testosterone pro...	1109024-02	60	1503607
Euro OTC	Tetracaine base	1512022	45	20160128*
Euro OTC	Tetracaine base	1405025-01	60	AR-14-FG-012475-01
Euro OTC	Tetracaine hydro...	1512019	45	20151228*
Euro OTC	Tetracaine hydro...	1807017	40	20180730*
Euro OTC	Tetracaine hydro...	1612005	40	20170112*
Euro OTC	Tetracaine hydro...	1506040	60	20150713*
Euro OTC	Tetracycline hyd...	1302012-01	40	1411585
Euro OTC	Tetracycline hyd...	1506012	60	20150707*
Euro OTC	Thiamine chlorid...	1211008-02	60	not necessary
Euro OTC	Thiamine chlorid...	1612004	40	20170117*
Euro OTC	Thymol	1312022-02	60	not necessary
Euro OTC	Thymol	1510004	60	20151016*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Tobramycin sulfate	1301011-02	60	1508071
Euro OTC	Tobramycin sulfate	1510035	60	20151208*
Euro OTC	Tobramycin sulfate	1510035-01	60	9202
Euro OTC	Tobramycin sulfate	1505015-01	60	8963
Euro OTC	Tobramycin sulfate	1807045	40	20180917*
Euro OTC	Tretinoin	1604059	40	20160603*
Euro OTC	Tretinoin	1402007-02	60	1507085
Euro OTC	Tretinoin	1402007-02	60	1507085
Euro OTC	Tretinoin	1510007	60	20151102*
Euro OTC	Tretinoin	1510007	60	20151102*
Euro OTC	Triamcinolone ac...	1608005	40	20160906*
Euro OTC	Triamcinolone ac...	1608005	22	20160906*
Euro OTC	Triamcinolone ac...	1601021	45	20160119*
Euro OTC	Triamcinolone ac...	1508008	60	20150902*
Euro OTC	Triamcinolone ac...	1605010	40	20160531*
Euro OTC	Triamcinolone ac...	1303001-02	40	not necessary
Euro OTC	Triamcinolone ac...	1605010	25	20160531*
Euro OTC	Triamcinolone ac...	1211013-02	30	not necessary
Euro OTC	Triclosan	1510045	60	20150511*
Euro OTC	Triclosan	1406032-01	60	not necessary
Euro OTC	Triclosan	1701032	40	20170207*
Euro OTC	Triclosan	1507041	60	20150804*
Euro OTC	Triclosan	1602003	45	20160204*
Euro OTC	Trometamol	1803016	40	20180320*
Euro OTC	Ubidecarenone	1703026	40	20170406*
Euro OTC	Ubidecarenone	1205029-01	40	not necessary
Euro OTC	Ubidecarenone	1605030	40	20160623*
Euro OTC	Ubidecarenone	1801019	40	20180126*
Euro OTC	Urea	1601044	40	20160208*
Euro OTC	Urea	1604044	40	20160512*
Euro OTC	Urea	1404024-01	60	not necessary
Euro OTC	Valine	1701026	40	20170222*
Euro OTC	Valine	1904006	40	20190423*
Euro OTC	Valine	1604055	40	20160601*
Euro OTC	Vancomycin hydro...	1908025	40	20190923*
Euro OTC	Vancomycin hydro...	1506029	60	AR-15-FG-009077-01
Euro OTC	Vancomycin hydro...	1610038	40	20161208*
Euro OTC	Vancomycin hydro...	1810013	40	20181029*
Euro OTC	Vancomycin hydro...	1407047-01	60	AR-15-FG-008950-01
Euro OTC	Zinc gluconate	1603035	40	20160427*
Euro OTC	Zinc gluconate	1807019	40	20180725*
Euro OTC	Zinc gluconate	1508007	60	20150827*
Euro OTC	Zinc gluconate	1407033-02	60	not necessary
Fagron	5-aminolevulinic...	17A12-B07-332960	40	not necessary
Fagron	5-aminolevulinic...	17A12-B07-332960	40	not necessary
Fagron	5-aminolevulinic...	13K27-N02	40	not necessary
Fagron	Acetylcysteine	15C16-B06-308039	40	not necessary
Fagron	Acetylsalicylic ...	14E02-B04-297099	40	not necessary
Fagron	Aciclovir	14G07-B05-296462	60	not necessary
Fagron	Aciclovir	13B04-B40-302333	60	not necessary
Fagron	Aciclovir	13B04-B40-294559	60	AR-15-FG-002570-01
Fagron	Albendazole	14B05-B41-299704	60	not necessary
Fagron	Albendazole	16B04-B06-322057	40	not necessary
Fagron	Albendazole	12J25-B40-293217	45	AR-15-FG-002571-01
Fagron	Alfatradiol	15J27-B08-323877	40	not necessary
Fagron	Alfatradiol	15J27-B08-323877	40	not necessary
Fagron	Alfatradiol	13G19-N06	40	1402683
Fagron	Alfatradiol	15J27-B08-323877	40	not necessary
Fagron	Allantoin	13I03-B05	40	not necessary
Fagron	Allopurinol	14E15-B05-307746	60	AR-16-FG-004367-01

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Allopurinol	14E15-B05-307746	40	not necessary
Fagron	Allopurinol	14E15-B05-315557	60	not necessary
Fagron	alpha-tocopheryl...	14B03-B04-295113	60	not necessary
Fagron	alpha-tocopheryl...	13H30-B02-287169	60	AR-15-FG-006388-01
Fagron	Ambroxol hydroch...	19E01-B02-361735	40	not necessary
Fagron	Ambroxol hydroch...	13I24-N02	40	not necessary
Fagron	Amifampridine	16A22-B04-317264	60	not necessary
Fagron	Amifampridine	13F05-N01	41	not necessary
Fagron	Amitriptyline hy...	16C31-B01-325540	40	not necessary
Fagron	Amitriptyline hy...	18A16-B06-348643	40	1909242
Fagron	Amphotericin B	13I16-N02	40	1508271
Fagron	Amphotericin B	13I16-N02	60	not necessary
Fagron	Amphotericin B	16L06-B04-333234	40	not necessary
Fagron	Arginine	13J28-B02	40	AR-15-FG-008045-01
Fagron	Arginine	16H24-B04-335673	40	not necessary
Fagron	Ascorbic acid	13G29-B07	40	1410500
Fagron	Ascorbyl palmitate	15F01-B99-309774	60	not necessary
Fagron	Ascorbyl palmitate	13K04-N01	40	1403226
Fagron	Ascorbyl palmitate	19D25-B04-361254	40	not necessary
Fagron	Ascorbyl palmitate	17E18-B06-335527	40	not necessary
Fagron	Ascorbyl palmitate	17E18-B06-335527	40	not necessary
Fagron	Ascorbyl palmitate	13C06-B40-291206	60	not necessary
Fagron	Atropine sulfate	16G20-B06-328437	40	not necessary
Fagron	Bacitracin	16B17-B03-322939	40	not necessary
Fagron	Bacitracin	15C12-B06-309185	60	9466
Fagron	Bacitracin	14E23-B03-296243	60	1503715
Fagron	Bacitracin	13A-03-N11	40	1401189
Fagron	Bacitracin	15C12-B06-309185	60	not necessary
Fagron	Bacitracin	14A21-B09-292909	60	1511741
Fagron	Beclometasone di...	13E02-B40-296353	29	not necessary
Fagron	Beclometasone di...	13I05-N03	60	1511689
Fagron	Beclometasone di...	13E02-B40-293328	40	not necessary
Fagron	Beclometasone di...	13E02-B40-296353	60	not necessary
Fagron	Benzocaine	16D13-B01-342411	40	1807429
Fagron	Benzocaine	14G02-B01-306980	40	not necessary
Fagron	Benzydamine hydr...	14I30-B04-300713	45	AR-15-FG-015977-01
Fagron	Benzydamine hydr...	15I17-B02-315060	60	not necessary
Fagron	Betaine monohydr...	14I26-B06-309878	60	not necessary
Fagron	Betaine monohydr...	18H13-B01-355585	40	not necessary
Fagron	Betamethasone di...	13I17-N05	80	not necessary
Fagron	Betamethasone di...	16C04-B03-324513	39	not necessary
Fagron	Betamethasone di...	13I17-N05	40	1401377
Fagron	Betamethasone di...	16C04-B03-324513	40	not necessary
Fagron	Betamethasone va...	12K07-N23	40	not necessary
Fagron	Betamethasone va...	13K29-B32-288243	60	not necessary
Fagron	Betamethasone va...	13K05-B08-292651	60	1511154
Fagron	Betamethasone, m...	12D26-N08	40	1402096
Fagron	Betamethasone, m...	11H31-B06-288219	40	1608324
Fagron	Betamethasone, m...	13H08-N06	60	1411037
Fagron	Bifonazole	16D07-B10-324381	40	not necessary
Fagron	Bifonazole	12G11-N09	40	not necessary
Fagron	Biotin	17K20-B03-351524	40	not necessary
Fagron	Biotin	13K29-N02	40	not necessary
Fagron	Bismuth gallath...	15B02-B04-309481	40	not necessary
Fagron	Bismuth gallath...	15J26-B02-317883	40	not necessary
Fagron	Boric acid	15L09-B07-322988	40	not necessary
Fagron	Boric acid	13I04-N06	40	not necessary
Fagron	Budesonide, micr...	16J04-B01-326581	40	not necessary
Fagron	Budesonide, micr...	16J04-B01-326581	30	not necessary
Fagron	Budesonide, micr...	18C15-B03-355399	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Budesonide, micr...	16J04-B01-326581	40	not necessary
Fagron	Calcium glycerop...	13I17-N04	40	not necessary
Fagron	Calcium glycerop...	15B03-B03-305050	60	1706313
Fagron	Calcium lactate ...	12E10-N03	40	not necessary
Fagron	Calcium lactate ...	12I12-N06	40	AR-15-FG-007939-01
Fagron	Calcium orotate ...	12L04-N07	40	not necessary
Fagron	Calcium orotate ...	13J23-N05	40	not necessary
Fagron	Calcium pantothe...	13H28-N11	40	not necessary
Fagron	Camphor, racemic	16E23-B08-336937	40	not necessary
Fagron	Camphor, racemic	14E19-B07-300365	60	not necessary
Fagron	Capsaicin	13J07-N07	40	AR-15-FG-007151-01
Fagron	Captopril	12G25-N08	40	not necessary
Fagron	Captopril	14B26-B01-312221	60	not necessary
Fagron	Carmellose sodiu...	18E15-B04-354617	40	not necessary
Fagron	Carmellose sodiu...	13D22-N02	40	1401187
Fagron	Cellulose, micro...	16B24-B11-323324	40	not necessary
Fagron	Cellulose, micro...	12C23-N06	70	1403080
Fagron	Cetomacrogol emu...	13C05-B40-289262	40	1408003
Fagron	Cetomacrogol emu...	12H07-N06	40	not necessary
Fagron	Cetostearyl alco...	13C13-N13	40	1402265
Fagron	Chloral hydrate	17E05-B05-335580	40	not necessary
Fagron	Chloral hydrate	13F06-N07	40	not necessary
Fagron	Chloramphenicol	13A11-N02	40	not necessary
Fagron	Chloramphenicol	14E01-B02-308559	60	not necessary
Fagron	Chloramphenicol	11F20-N05	120	1412249
Fagron	Chlorhexidine di...	15J02-B08-322916	40	not necessary
Fagron	Chlorhexidine di...	14I08-B02-306223	60	not necessary
Fagron	Chlorhexidine di...	14I08-B02-306223	40	not necessary
Fagron	Chlorhexidine di...	15J02-B08-322916	40	not necessary
Fagron	Chlortetracyclin...	12K30-N07	40	1405318
Fagron	Chlortetracyclin...	14J08-B07-307742	60	not necessary
Fagron	Chlortetracyclin...	12K30-N07	40	not necessary
Fagron	Cholesterol	17E17-B05-339074	40	not necessary
Fagron	Cholesterol	13F10-B40	40	not necessary
Fagron	Cholesterol	15J23-B03-325106	40	not necessary
Fagron	Ciclopirox olamine	14F13-B03-302340	60	1608170
Fagron	Ciclopirox olamine	16G26-B02-324884	40	not necessary
Fagron	Ciclopirox olamine	14F13-B03-302340	60	1608170
Fagron	Ciclopirox olamine	13A16-N04	40	not necessary
Fagron	Ciclopirox olamine	16G26-B02-324884	40	not necessary
Fagron	Ciclosporin	13D26-N07	42	1508429
Fagron	Ciclosporin	18G30_B04-354759	40	not necessary
Fagron	Ciclosporin	13D26-N07	40	1402057
Fagron	Ciprofloxacin hy...	12B08-N03	40	1401421
Fagron	Ciprofloxacin hy...	14B25-B02	40	not necessary
Fagron	Ciprofloxacin hy...	18F27-B01-353346	40	not necessary
Fagron	Citric acid	13G24-N13	40	not necessary
Fagron	Citric acid	12F18-N03	40	not necessary
Fagron	Citric acid mono...	12D26-N02	40	not necessary
Fagron	Citrulline	18K15-B01-356962	40	not necessary
Fagron	Citrulline	16D11-B02-326797	40	not necessary
Fagron	Clindamycin hydr...	16D05-B08-332106	40	not necessary
Fagron	Clindamycin hydr...	13G24-N09	40	1412075
Fagron	Clioquinol	14I11-B03-300049	60	not necessary
Fagron	Clioquinol	14I11-B03-300049	60	not necessary
Fagron	Clioquinol	13I26-N15	40	not necessary
Fagron	Clobetasol propi...	18I12-B02-357669	40	not necessary
Fagron	Clonidine hydroc...	14H11-B04-311163	60	AR-18-FG-006797-02
Fagron	Clonidine hydroc...	19B27-B02-361777	40	not necessary
Fagron	Clonidine hydroc...	14H11-B04-311163	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Clonidine hydroc...	13I25-B07	40	not necessary
Fagron	Clonidine hydroc...	14H11-B04-311163	40	not necessary
Fagron	Clotrimazole	14C04-B31	40	not necessary
Fagron	Clotrimazole	14B05-B36	40	not necessary
Fagron	Clotrimazole	12I21-N01	40	not necessary
Fagron	Colistin sulfate	14E02-B03	40	1501654
Fagron	Colistin sulfate	13G05-N01	40	not necessary
Fagron	Colistin sulfate	16F03-B01-324216	40	not necessary
Fagron	Colistin sulfate	13G05-N01	40	not necessary
Fagron	Cyproterone acet...	13J16-No4	75	1602654
Fagron	Cyproterone acet...	13C11-N02	40	1501534
Fagron	Cyproterone acet...	17D11-B07-333693	20	not necessary
Fagron	Cyproterone acet...	17D11-B07-333693	40	not necessary
Fagron	Cyproterone acet...	15I11-B01-322117	40	not necessary
Fagron	Cystine	15F08-B05-310685	40	not necessary
Fagron	Deoxycholic acid	14A02-B32-303429	60	AR-16-FG-004648-01
Fagron	Dequalinium chlo...	15B16-B04-332113	40	not necessary
Fagron	Dequalinium chlo...	13I26-N10	40	not necessary
Fagron	Dexamethasone	16F02-B05-324844	40	not necessary
Fagron	Dexamethasone	19A24-B01-361445	40	not necessary
Fagron	Dextromethorphan...	15D23-B07-308849	60	AR-16-FG-008487-02
Fagron	Dextromethorphan...	16A22-B03-321371	40	not necessary
Fagron	Dextromethorphan...	17B28-B01-341594	40	not necessary
Fagron	Dextromethorphan...	15D23-B07-308849	60	not necessary
Fagron	D-galactose	14H05-B09-310972	60	not necessary
Fagron	D-galactose	13I13-B40-287993	60	AR-15-FG-008936-01
Fagron	Diclofenac sodium	13E03-N02	40	not necessary
Fagron	Diclofenac sodium	12L28-B05	40	not necessary
Fagron	Diltiazem hydroc...	16G20-B03-328133	40	not necessary
Fagron	Diltiazem hydroc...	16G20-B03-328133	40	not necessary
Fagron	Diltiazem hydroc...	14C12-B04	40	not necessary
Fagron	Diphenhydramine ...	18K28-B02-358954	40	not necessary
Fagron	Diphenhydramine ...	16C25-B03-331026	40	AR-17-FG-015432-02
Fagron	Diphenylcyclopro...	18I03-B04-356833	40	not necessary
Fagron	Diphenylcyclopro...	15J12-B01-316620	60	9341
Fagron	Diphenylcyclopro...	12L17-N07	40	1407599
Fagron	Diphenylcyclopro...	16F02-B10-329708	40	1902476
Fagron	Dithranol	18B23-B05-348137	40	not necessary
Fagron	Doxycycline hycl...	12G12-N15	40	not necessary
Fagron	Doxycycline hycl...	12G12-N15	40	1404233
Fagron	Doxycycline hycl...	14K28-B05-312218	60	not necessary
Fagron	Doxylamine hydro...	15F01-B04-309189	40	AR-16-FG-011361-01
Fagron	Doxylamine hydro...	14F03-B01-293498	60	1410524
Fagron	Emulsifying ceto...	16I13-B01-326511	40	not necessary
Fagron	Emulsifying ceto...	13L17-B05	40	1505070
Fagron	Eosin Y	13F14-T01-017395	40	1608469
Fagron	Eosin Y	13F14-T01-017395	40	not necessary
Fagron	Eosin Y	17K21-T07-042542	40	1805338
Fagron	Eosin Y	13F14-T01-017395	40	1801115
Fagron	Eosin Y	12E25-N01	60	not necessary
Fagron	Eosin Y	13F14-T01-017395	40	1705197
Fagron	Eosin Y	15G01-T03-025958	60	1604114
Fagron	Ephedrine hydroc...	15H28-B02-316349	40	not necessary
Fagron	Ephedrine hydroc...	17G05-B05-337987	40	not necessary
Fagron	Epinephrine hydr...	12L22-N04	40	not necessary
Fagron	Epinephrine hydr...	18K30-B03-358514	40	not necessary
Fagron	Ergotamine tartr...	15B16-B05-308418	45	9384
Fagron	Ergotamine tartr...	15B16-B05-308417	40	not necessary
Fagron	Ergotamine tartr...	15B16-B05-308418	45	9384
Fagron	Ergotamine tartr...	15B16-B05-308417	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Erythromycin	12E24-N04	60	not necessary
Fagron	Erythromycin	14B06-B35-290872	60	not necessary
Fagron	Escin soluble	12I07-B40-306607	40	1610369
Fagron	Escin soluble	12K22-N13	40	not necessary
Fagron	Estradiol benzoate	12H02-N02	40	not necessary
Fagron	Estradiol benzoate	14B20-B02	40	not necessary
Fagron	Estradiol benzoate	14B20-B02	60	not necessary
Fagron	Estradiol benzoate	12H02-N02	60	not necessary
Fagron	Estradiol hemihy...	14G24-B04-297404	55	not necessary
Fagron	Estradiol hemihy...	16H04-B01-328738	40	not necessary
Fagron	Estradiol valerate	12D18-N02	40	not necessary
Fagron	Estradiol valerate	14E26-B02-294846	60	not necessary
Fagron	Estradiol valerate	12D18-N02	40	not necessary
Fagron	Estriol	16D15-B03-322957	40	not necessary
Fagron	Estriol	16D15-B03-322957	40	not necessary
Fagron	Estriol	11J13-N07	40	1412180
Fagron	Estriol	16I19-B01-331468	40	not necessary
Fagron	Ethacridine lact...	15A29-B02-306219	40	not necessary
Fagron	Ethacridine lact...	13I24-N03	40	AR-15-FG-007011-02
Fagron	Ethinylestradiol	14E15-B03-296049	55	1507109
Fagron	Ethinylestradiol	14E15-B03-296049	56	1507109
Fagron	Ethinylestradiol	14E15-B03-310270	60	not necessary
Fagron	Ethinylestradiol	13I09-N04	40	1402058
Fagron	Ethinylestradiol	14E15-B03-310270	60	not necessary
Fagron	Ethinylestradiol	14E15-B03-310270	60	not necessary
Fagron	Ethinylestradiol	18E25-B01-355721	40	not necessary
Fagron	Ethyl hydrogen f...	13B25-N08	60	1411090
Fagron	Ethyl hydrogen f...	13G10-N01	60	1411091
Fagron	Fampridine	13G26-N06	40	1402686
Fagron	Fampridine	17K20-B02-349165	40	not necessary
Fagron	Fampridine	16D11-B07-318220	40	AR-16-FG-011566-01
Fagron	Fampridine	14D11-B06-293750	60	not necessary
Fagron	Fampridine	16E25-B03-324097	40	not necessary
Fagron	Fampridine	17K20-B02-349165	40	not necessary
Fagron	Fampridine	16E25-B03-324097	40	not necessary
Fagron	Finasteride (for...	13I03-N01	40	1411272
Fagron	Finasteride (for...	14A16-B40-295980	45	not necessary
Fagron	Finasteride (for...	14A16-B40-295980	60	not necessary
Fagron	Fluocinolonaceto...	14L19-B04-307101	40	not necessary
Fagron	Fluocinolonaceto...	14L19-B04-307101	45	not necessary
Fagron	Fluocinolonaceto...	17F14-B08-350105	40	not necessary
Fagron	Fluocinolonaceto...	14L19-B04-307101	45	not necessary
Fagron	Folic acid hydrate	12G17-N05	40	not necessary
Fagron	Folic acid hydrate	16C03-B06-322513	40	not necessary
Fagron	Fructose	16H05-B01-326259	40	not necessary
Fagron	Fructose	16H05-B01-326259	40	not necessary
Fagron	Furosemide	13A22-N08	40	1401379
Fagron	Furosemide	17D06-B01-338982	40	not necessary
Fagron	Furosemide	14K10-B02-301746	60	not necessary
Fagron	Fusidic acid	13L13-B01-289513	60	not necessary
Fagron	Fusidic acid	14E01-B08-299791	60	not necessary
Fagron	Fusidic acid	12G04-N01	40	1401162
Fagron	Gentamicin sulfate	16B29-B01-328731	40	not necessary
Fagron	Gentamicin sulfate	14E08-B03-294848	60	not necessary
Fagron	Gentamicin sulfate	12B21-N02	40	1411075
Fagron	Glucose monohydr...	14B04-B01	40	not necessary
Fagron	Glucose monohydr...	11L19-N13	40	not necessary
Fagron	Glutamic acid	15G31-B08-311884	40	not necessary
Fagron	Glutamic acid	11I30-N02	40	AR-15-FG-010446-02
Fagron	Glutamine	17E05-B06-334869	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Glutamine	17E05-B06-334869	40	AR-17-FG-016178-01
Fagron	Glycerol monoste...	15E13-B03-305532	60	not necessary
Fagron	Glycerol monoste...	14C12-B40-292622	60	not necessary
Fagron	Glycine	13G09-N03	40	not necessary
Fagron	Glycine	12G13-N04	40	1401163
Fagron	Glycopyrronium b...	18G25-B01-356318	40	not necessary
Fagron	Glycopyrronium b...	16E09-B08-331027	40	1707615
Fagron	Glycopyrronium b...	16E09-B08-331027	40	not necessary
Fagron	Glycopyrronium b...	18C30-B07-350744	40	not necessary
Fagron	Gramicidin	18A19-B14-350428	40	not necessary
Fagron	Gramicidin	18L20-B01-358528	40	not necessary
Fagron	Gramicidin	13J15-N02	40	AR-15-FG-000813-01
Fagron	Gramicidin	17E09-B02-341211	40	AR-18-FG-008970-02
Fagron	Griseofulvin	13D19-N03	40	1401419
Fagron	Griseofulvin	16F13-B06-330832	40	not necessary
Fagron	Griseofulvin	14E16-B01-296834	60	not necessary
Fagron	Hard fat W35	15A07-B12-302751	60	1604639
Fagron	Hard fat W35	15A07-B12-309245	60	not necessary
Fagron	Hard fat W35	16K08-B10-334910	40	not necessary
Fagron	Heparin sodium	12H30-N07	40	1404230
Fagron	Heparin sodium	12L20-B01-317990	40	not necessary
Fagron	Histamine dihydr...	19D19-B02-361228	40	not necessary
Fagron	Histamine dihydr...	14G16-B02-299175	100	not necessary
Fagron	Histidine	18J03-B02-355093	40	not necessary
Fagron	Histidine	15I01-B12-329082	40	not necessary
Fagron	Hydrochlorothiaz...	17L20-B02-345079	40	not necessary
Fagron	Hydrochlorothiaz...	11L01-N02	40	not necessary
Fagron	Hydrocortisone	16H01-B01-326071	40	not necessary
Fagron	Hydrocortisone	16A19-B06-321981	40	1702302
Fagron	Hydrocortisone	16A19-B06-320970	30	not necessary
Fagron	Hydrocortisone	16A19-B06-321981	40	not necessary
Fagron	Hydrocortisone	16H01-B01-326071	40	not necessary
Fagron	Hydrocortisone	16A19-B06-320970	40	not necessary
Fagron	Hydrocortisone a...	18H07-B03-356259	40	not necessary
Fagron	Hydrocortisone a...	15A30-B09-307090	60	not necessary
Fagron	Hydrocortisone b...	15J12-B03-322394	40	not necessary
Fagron	Hydroxyethylcell...	13H22-N11	40	not necessary
Fagron	Hydroxyethylcell...	16E09-B05-327042	40	not necessary
Fagron	Hyoscine butylbr...	13K22-N03	40	not necessary
Fagron	Hyoscine butylbr...	12K01-N07	40	not necessary
Fagron	Hyoscine butylbr...	15I09-B04-312245	60	not necessary
Fagron	Hyoscine butylbr...	12K01-N07	40	1406057
Fagron	Hyoscine butylbr...	18A03-B01-352889	40	not necessary
Fagron	Hyoscine hydrobr...	15L09-B04-334706	40	not necessary
Fagron	Hyoscine hydrobr...	15L09-B04-334706	40	not necessary
Fagron	Hypromellose 4000	15G08-B14-314845	40	not necessary
Fagron	Inositol	13D15-B40-328250	40	AR-17-FG-016393-01
Fagron	Inositol	13D15-B40-328250	40	not necessary
Fagron	Iodoform	12L11-B41	40	not necessary
Fagron	Ipratropium brom...	12K01-N08	40	not necessary
Fagron	Ipratropium brom...	12C15-N01	40	1403675
Fagron	Isosorbide dinit...	15I10-B01-328577	30	not necessary
Fagron	Isosorbide dinit...	15I10-B01-328577	40	1709092
Fagron	Isosorbide dinit...	17C29-B01-339447	40	not necessary
Fagron	Isosorbide dinit...	13D15-N02	39	not necessary
Fagron	Isosorbide dinit...	12F08-N04	40	not necessary
Fagron	Ketamine hydroch...	13H19-N01	40	8931
Fagron	Ketamine hydroch...	14I03-B02-303165	60	not necessary
Fagron	Ketamine hydroch...	14I03-B02-311706	40	not necessary
Fagron	Ketoconazole	16A26-B01-322982	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Ketoconazole	12A16-N06	40	not necessary
Fagron	Levocarnitine	12H10-N07	40	not necessary
Fagron	Levocarnitine	14A20-B11-290728	60	AR-15-FG-016322-01
Fagron	Levocarnitine	14A20-B11-290728	60	AR-16-FG-005969-01
Fagron	Levocarnitine	16L06-B03-338285	40	not necessary
Fagron	Levocarnitine	15L04-B05-319229	40	AR-17-FG-012937-01
Fagron	Levothyroxin sod...	14E28-B05-306102	60	not necessary
Fagron	Levothyroxin sod...	14E28-B05-306102	53	not necessary
Fagron	Levothyroxin sod...	14E28-B05-298322	47	1506426
Fagron	Levothyroxin sod...	18G27-B02-352032	40	not necessary
Fagron	Levothyroxin sod...	14E28-B05-298322	30	1506426
Fagron	L-glutathione	18C01-B08-347604	40	not necessary
Fagron	Lidocaine	13H05-N05	40	1401378
Fagron	Lidocaine	16F21-B06-327232	40	not necessary
Fagron	L-isoleucine	13E02-B01-286930	40	1403154
Fagron	L-isoleucine	13E02-B01	40	not necessary
Fagron	Lithium orotate ...	15F25-B01-315219	40	not necessary
Fagron	Lithium orotate ...	13J04-N01	40	1405616
Fagron	Lithium orotate ...	13I09-B04-290311	60	not necessary
Fagron	L-leucine	13J09-N03	40	1402431
Fagron	L-leucine	15D09-B06-322249	40	not necessary
Fagron	L-malic acid	15J09-B03-313987	40	not necessary
Fagron	L-malic acid	15J09-B03-313987	40	not necessary
Fagron	L-malic acid	16J06-B98-327269	40	not necessary
Fagron	L-malic acid	15J09-B03-313987	60	9346
Fagron	L-malic acid	16J06-B98-327269	40	not necessary
Fagron	L-phenylalanine	17B13-B03-343288	40	not necessary
Fagron	L-proline	16J14-B01-327195	40	not necessary
Fagron	L-proline	15D13-B02-308152	60	1704059
Fagron	L-threonine	12H02-N03	40	1402580
Fagron	L-tryptophan	17G19-B01-343931	40	not necessary
Fagron	L-tryptophan	14B05-B32-293485	60	AR-15-FG-009473-01
Fagron	Lysine hydrochloride...	14H20-B03-301894	60	AR-16-FG-005856-01
Fagron	Lysine hydrochloride...	15B09-B06-304170	60	not necessary
Fagron	Macrogols 4000	16F10-B04-325645	40	not necessary
Fagron	Macrogols 4000	15B13-B04-305461	60	AR-16-FG-004020-01
Fagron	Magnesium stearate	12C29-N06	110	not necessary
Fagron	Magnesium stearate	12D26-B03	40	1407310
Fagron	Mannitol	15A26-B01-305841	40	not necessary
Fagron	Mannitol	17B08-B04-352402	40	not necessary
Fagron	Mannitol	15A26-B01-305841	40	not necessary
Fagron	Melatonin	15G06-B03-323733	40	not necessary
Fagron	Melatonin	15G06-B03-315183	40	not necessary
Fagron	Menthol	18L04-T06-065445	40	not necessary
Fagron	Menthol	16L20-T04-044769	40	not necessary
Fagron	Mepivacaine hydr...	13L04-B07-297077	60	AR-15-FG-007876-01
Fagron	Mepivacaine hydr...	12A20-N12	40	1401555
Fagron	Mesalazine	15H13-B06-317177	60	not necessary
Fagron	Mesalazine	16I28-B02-327679	40	not necessary
Fagron	Methacholine chl...	12F01-N04	60	1403251
Fagron	Methacholine chl...	16L20-B01-329966	40	not necessary
Fagron	Methacholine chl...	16L20-B01-329968	40	not necessary
Fagron	Methionine	12K01-N03	40	1402432
Fagron	Methionine	13L19-B36-295115	60	not necessary
Fagron	Methoxsalen	18A12-B50-355102	40	not necessary
Fagron	Methoxsalen	15J15-B03-321321	40	not necessary
Fagron	Methyl-4-hydroxy...	14L30-T01-018345	40	not necessary
Fagron	Methyl-4-hydroxy...	19B05-T03-062158	40	not necessary
Fagron	Methylcellulose	12D02-N09	40	1402447
Fagron	Methylcellulose	15H04-B07-331291	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Methylprednisolone	18E16-B02-B02-358166	40	not necessary
Fagron	Methylprednisolone	15I04-B05-311572	45	not necessary
Fagron	Methylprednisolone	15I04-B05-311572	45	not necessary
Fagron	Methylprednisolone	15I04-B05-311572	45	AR-16-FG-009102-01
Fagron	Methylprednisolone	15I04-B05-311572	45	AR-16-FG-009102-01
Fagron	Methylprednisolone	13G19-N05	40	1402438
Fagron	Metoprolol tartr...	16E31-B03-327107	40	not necessary
Fagron	Metoprolol tartr...	16E31-B03-327107	40	not necessary
Fagron	Metoprolol tartr...	16E31-B03-327107	40	not necessary
Fagron	Metronidazole	16I19-B04-326155	40	not necessary
Fagron	Metronidazole	13D10-N02	40	not necessary
Fagron	Miconazole nitrate	18C23-B02-356141	40	not necessary
Fagron	Miconazole nitrate	13D22-N14	100	not necessary
Fagron	Mometasone furoa...	13F17-N01	44	not necessary
Fagron	Mometasone furoa...	13F17-N01	40	not necessary
Fagron	Mometasone furoa...	16A20-B06-319750	30	not necessary
Fagron	Mometasone furoa...	16A20-B06-319750	40	not necessary
Fagron	Mupirocin calcium	15H13-B01-316111	40	not necessary
Fagron	Mupirocin calcium	12E02-N07	40	1401556
Fagron	Naloxone hydroch...	13D22-N07	40	1402059
Fagron	Naloxone hydroch...	13L12-B03	40	not necessary
Fagron	Naphazoline nitr...	12L10-N03	40	not necessary
Fagron	Naphazoline nitr...	12J08-B40-287904	60	not necessary
Fagron	Naproxen	16B04-B09-326261	40	not necessary
Fagron	Naproxen	15F18-B03-315691	45	not necessary
Fagron	Naproxen	15F18-B03-309492	45	AR-16-FG-008486-02
Fagron	Neomycin sulfate	14C27-B05-302517	60	not necessary
Fagron	Neomycin sulfate	14C27-B05-302516	60	1604527
Fagron	Nicotinic acid	16H19-B02-330121	40	AR-17-FG-016013-01
Fagron	Nicotinic acid	14E01-B09-293947	60	AR-15-FG-008938-01
Fagron	Norethisterone a...	17B21-B01-341608	40	not necessary
Fagron	Norethisterone a...	12L14-N03	40	1401558
Fagron	Nystatin	19B28-B05-191210	40	not necessary
Fagron	Ofloxacin	15H19-B01-312926	60	not necessary
Fagron	Ofloxacin	15H19-B01-325226	40	1707611
Fagron	Ofloxacin	18J30-B01-360431	40	not necessary
Fagron	Ofloxacin	13A17-N01	60	not necessary
Fagron	Ofloxacin	13A17-N01	40	not necessary
Fagron	Ofloxacin	13A17-N01	40	1403155
Fagron	Omeprazole	15K12-B02-314608	40	not necessary
Fagron	Omeprazole	15K12-B02-314608	40	not necessary
Fagron	Omeprazole	19F14-B01-362735	40	1910227
Fagron	Ornithine aspart...	14I24-B04-307758	40	not necessary
Fagron	Ornithine aspart...	14I24-B04-307758	40	not necessary
Fagron	Ornithine aspart...	17I21-B05-342561	40	1804455
Fagron	Ornithine aspart...	14A20-B07-296093	45	1506279
Fagron	Ornithine hydroc...	13A03-N05	40	1403094
Fagron	Oxybuprocaine hy...	13D16-N12	40	1402060
Fagron	Oxybuprocaine hy...	15A26-B07-306787	60	not necessary
Fagron	Oxybuprocaine hy...	13J24-B07-295015	60	AR-15-FG-014812-02
Fagron	Oxybuprocaine hy...	17G21-B03-338576	40	not necessary
Fagron	Oxybuprocaine hy...	16L05-B06-334330	40	not necessary
Fagron	Oxybutynin hydro...	14I16-B03-312059	40	not necessary
Fagron	Oxybutynin hydro...	14A29-B09-293841	45	not necessary
Fagron	Oxybutynin hydro...	13H26-B03-290572	60	AR-15-FG-002173-01
Fagron	Oxytetracycline ...	15D29-B99-309066	60	not necessary
Fagron	Oxytetracycline ...	13D08-N11	60	not necessary
Fagron	Paracetamol	13G23-N03	40	1402297
Fagron	Paracetamol	13D16-N08	40	1402296
Fagron	Paracetamol	13D30-N03	40	1402271

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Phenazone	13G16-B07-287172	40	not necessary
Fagron	Phenazone	13G16-B07-287172	40	AR-16-FG-011360-01
Fagron	Phenazone	15H28-B09-322944	40	not necessary
Fagron	Pheniramine male...	13D24-N04	40	not necessary
Fagron	Phenylephrine hy...	16A20-B07-317987	40	AR-16-FG-013385-01
Fagron	Phenylephrine hy...	12J19-B40-288254	40	not necessary
Fagron	Phenylephrine hy...	14H13-B02-308438	60	AR-16-FG-008920-01
Fagron	Phenytoin	13G26-N02	40	1401380
Fagron	Phenytoin	18F01-B02-350492	40	not necessary
Fagron	Phenytoin	15E12-B05-312661	45	not necessary
Fagron	Phenytoin	18F01-B02-350492	40	not necessary
Fagron	Poloxamers 407	14G17-B03-296866	60	not necessary
Fagron	Polymyxin B sulf...	12G13-N05	40	1410347
Fagron	Polymyxin B sulf...	14F04-B03-296722	60	not necessary
Fagron	Polymyxin B sulf...	13I13-N05	45	AR-15-FG-001516-01
Fagron	Polymyxin B sulf...	13C12-N02	40	1402691
Fagron	Polymyxin B sulf...	15H14-B01-314364	60	not necessary
Fagron	Povidone 25, sol...	13J21-N02	40	AR-16-FG-019116-01
Fagron	Povidone 25, sol...	12E22-N05	40	1411537
Fagron	Povidone 25, sol...	15C12-B02-323006	40	1710348
Fagron	Prasterone	16L13-B01-330724	40	not necessary
Fagron	Prasterone	13J10-B02-287175	40	1404309
Fagron	Prasterone	15G14-B07-310827	60	not necessary
Fagron	Prednicarbate, m...	12D27-N04	48	not necessary
Fagron	Prednicarbate, m...	14A27-B02-291256	60	AR-16-FG-011442-01
Fagron	Prednicarbate, m...	12D27-N04	40	not necessary
Fagron	Prednisolon dihy...	18D23-B08-351348	40	1812421
Fagron	Prednisolon dihy...	14E07-B06-295251	60	not necessary
Fagron	Prednisolone ace...	12H14-N07	40	not necessary
Fagron	Prednisolone ace...	15I23-B04-345785	40	not necessary
Fagron	Prednisolone ace...	12H14-N07	49	not necessary
Fagron	Prednisolone hyd...	14G29-B07-304733	40	not necessary
Fagron	Prednisolone hyd...	14G29-B07-304733	40	1701370
Fagron	Prednisolone hyd...	13D16-N11	40	1509250
Fagron	Prednisolone, mi...	13H07-B02-285592	60	not necessary
Fagron	Prednisolone, mi...	13H07-B02-285592	60	not necessary
Fagron	Prednisolone, mi...	16C15-B01-320761	36	not necessary
Fagron	Prednisolone, mi...	14I10-B03-302263	40	not necessary
Fagron	Prednisolone, mi...	14I10-B03-302263	40	not necessary
Fagron	Prednisolone, mi...	16C15-B01-320761	40	not necessary
Fagron	Prednisone	13C14-B40-289148	60	not necessary
Fagron	Prednisone	16J19-B02-328642	40	not necessary
Fagron	Pregnenolone	15B11-B07-310278	45	9338
Fagron	Pregnenolone	18L27-B02-357840	40	not necessary
Fagron	Prilocaine hydro...	13I02-N03	40	not necessary
Fagron	Prilocaine hydro...	14H13-B03-297591	60	AR-15-FG-016323-01
Fagron	Prilocaine hydro...	13I02-N03	40	1403257
Fagron	Procaine hydroch...	15J27-B10-322937	40	not necessary
Fagron	Procaine hydroch...	13F12-B40-301798	60	not necessary
Fagron	Progesterone, mi...	16A20-B02-319103	40	not necessary
Fagron	Progesterone, mi...	16A20-B02-319103	40	not necessary
Fagron	Propranolol hydr...	13I13-N04	40	1402307
Fagron	Propranolol hydr...	15F16-B02-312979	40	not necessary
Fagron	Propyphenazone	16E02-B08-327928	40	not necessary
Fagron	Propyphenazone	16H25-B01-334972	40	not necessary
Fagron	Quinidine sulfat...	14I04-B09-304357	60	AR-15-FG-015409-01
Fagron	Quinidine sulfat...	12J16-B03	39	1412323
Fagron	Quinidine sulfat...	12L04-N10	40	AR-14-FG-013249-01
Fagron	Quinine hydrochl...	13I30-B02-295630	60	AR-15-FG-015410-01
Fagron	Quinine hydrochl...	13G10-N07	30	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Quinine sulfate ...	14B27-B01	40	AR-16-FG-009455-01
Fagron	Quinine sulfate ...	13J23-N03	40	1404401
Fagron	Resocrinol	12G24-N04	40	1401462
Fagron	Resocrinol	13I18-B09-291533	45	not necessary
Fagron	Riboflavin	16B18-B06-319026	40	AR-16-FG-014466-01
Fagron	Riboflavin	15G20-B13-310997	60	not necessary
Fagron	Riboflavin	17G14-B05-337885	40	AR-18-FG-000141-01
Fagron	Riboflavin	13D02-N03	40	not necessary
Fagron	Riboflavin	13D02-N03	40	AR-15-FG-000733-01
Fagron	Riboflavin phosph...	17A03-B01-335897	40	AR-18-FG-000920-01
Fagron	Riboflavin phosph...	13I26-N06	40	1403705
Fagron	Riboflavin phosph...	15J02-B03-313137	60	AR-16-FG-009883-02
Fagron	Rice starch	16B25-B01-327231	40	not necessary
Fagron	Rice starch	16B25-B01-327231	40	AR-17-FG-015298-02
Fagron	Saccharin sodium	13I20-B06-291654	40	AR-17-FG-012457-01
Fagron	Saccharin sodium	14F30-B02-296322	40	AR-17-FG-022421-01
Fagron	Saccharin sodium	14F30-B02-296322	40	AR-18-FG-007117-02
Fagron	Salicylic acid	16I15-B08-333932	40	not necessary
Fagron	Salicylic acid	190263	40	not necessary
Fagron	Sildenafil citrate	13J15-N01	80	not necessary
Fagron	Sildenafil citrate	13F24-N03	40	1402052
Fagron	Sodium citrate	14A13-B02-294541	60	not necessary
Fagron	Sodium citrate	12E04-N05	40	not necessary
Fagron	Sodium citrate	16F10-B01-326088	40	not necessary
Fagron	Sodium cromoglic...	12C09-N02	79	1401553
Fagron	Sodium cromoglic...	13G24-N04	40	not necessary
Fagron	Sodium edetate	16F21-B02-325085	40	not necessary
Fagron	Sodium edetate	13E15-N05	40	1402193
Fagron	Sorbitan stearate	14J07-B02-299399	60	AR-16-FG-004368-01
Fagron	Spironolactone	2294216E20-B06-322016	40	not necessary
Fagron	Spironolactone	15G22-B10-311159	40	not necessary
Fagron	SyrSpend® SF Alka	18C13-B11-350851	40	not necessary
Fagron	SyrSpend® SF Alka	17D11-B04-332593	40	not necessary
Fagron	SyrSpend® SF Alka	17F14-B02-335945	20	1804106
Fagron	SyrSpend® SF Alka	17D11-B04-332593	40	not necessary
Fagron	SyrSpend® SF Alka	18B21-B08-347192	40	1809264
Fagron	SyrSpend® SF Alka	15H11-B04-311722	40	not necessary
Fagron	SyrSpend® SF pH4	15A26-B04-305667	40	1701055
Fagron	SyrSpend® SF pH4	15A26-B04-311720	40	not necessary
Fagron	SyrSpend® SF pH4	15A26-B04-311720	40	1608352
Fagron	Tannin	13H30-B03-287645	60	1711482
Fagron	Tannin	16D08-B04-328410	40	AR-17-FG-015963-01
Fagron	Tannin	16D08-B04-328410	40	not necessary
Fagron	Tartaric acid	13I25-B01-294093	60	1711509
Fagron	Tartaric acid	14G07-B02-308554	40	not necessary
Fagron	Testosterone	17J26-B01-342032	40	not necessary
Fagron	Testosterone	12F13-N02	40	not necessary
Fagron	Testosterone	15E18-B08-307772	54	not necessary
Fagron	Testosterone	15E18-B08-307772	60	not necessary
Fagron	Tetracaine base	14A16-B04-292916	60	not necessary
Fagron	Tetracaine base	12L19-N10	40	not necessary
Fagron	Tetracaine base	14A16-B04-292915	60	not necessary
Fagron	Tetracaine hydro...	16G28-B03-327368	40	not necessary
Fagron	Tetracaine hydro...	15K04-B03-322151	40	not necessary
Fagron	Thalidomide	15A13-B01-312230	40	1704509
Fagron	Thalidomide	14E27B04295348	60	1503608
Fagron	Thymol	15H13-B05-328912	40	not necessary
Fagron	Tobramycin sulfate	15I17-B04-315697	45	not necessary
Fagron	Tobramycin sulfate	12G03-N03	40	1404541
Fagron	Tobramycin sulfate	16B05-B03-336397	40	1801385

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Tranexamic acid	14L15-B06-3000036	60	AR-15-FG-008937-02
Fagron	Tranexamic acid	14C21-B08-292225	60	AR-17-FG-016394-01
Fagron	Tretinoin	15F01-B06-317884	40	not necessary
Fagron	Tretinoin	15F01-B06-317884	40	not necessary
Fagron	Triamcinolone ac...	12E15-N03	60	1412051
Fagron	Triamcinolone ac...	15C19-B06-308572	60	not necessary
Fagron	Triamcinolone, m...	18A19-B05-351174	32	1811263
Fagron	Triamcinolone, m...	18E09-B04-353921	40	1907730
Fagron	Triamcinolone, m...	12A03-N09	40	not necessary
Fagron	Triamcinolone, m...	18A19-B05-351174	40	1811263
Fagron	Triamcinolone, m...	14B13-B30-316101	40	not necessary
Fagron	Triamcinolone, m...	12K07-N13	80	AR-15-FG-003373-01
Fagron	Triclosan	17K02-B06-343548	40	not necessary
Fagron	Trometamol	13E22-N05	40	AR-17-FG-016395-01
Fagron	Trometamol	18A23-B10-351258	40	not necessary
Fagron	Tropicamide	14D23-B04-293294	60	AR-15-FG-001641-01
Fagron	Tropicamide	16G19-B05-324829	40	not necessary
Fagron	Tropicamide	17J06-B04-341322	40	1807428
Fagron	Tyrosine	18E15-B02-359934	40	not necessary
Fagron	Tyrosine	13G24-N22	60	1411030
Fagron	Tyrosine	16B26-B01-322911	40	not necessary
Fagron	Tyrothricin	12K14-N01	40	1509178
Fagron	Tyrothricin	12K14-N01	60	AR-15-FG-001172-01
Fagron	Tyrothricin	14A23-B05-294146	60	AR-15-FG-015005-01
Fagron	Ubidecarenone	16G26-B01-331742	40	not necessary
Fagron	Ubidecarenone	18J24-B01-357225	40	not necessary
Fagron	Urea	14I18-B03-303581	60	not necessary
Fagron	Urea	12F25-N05	40	not necessary
Fagron	Urea	14I18-B03-303581	60	not necessary
Fagron	Uridine	11G07-N04	60	1402481
Fagron	Uridine	13K18-B31-288497	60	not necessary
Fagron	Valine	15A14-B07-307756	60	9198
Fagron	Vancomycin hydro...	18L18-B02-359948	40	not necessary
Fagron	Vancomycin hydro...	17I18-B08-340944	40	not necessary
Fagron	Vanillin	13J01-N07	60	not necessary
Fagron	Xylitol	14G04-B04-298356	60	AR-15-FG-008046-01
Fagron	Xylitol	14G04-B04-303524	60	not necessary
Fagron	Xylometazoline h...	14K19-B05-305006	60	not necessary
Fagron	Xylometazoline h...	13I18-B03-295481	60	AR-15-FG-007785-01
Fagron	Xylometazoline h...	14K19-B05-305006	60	not necessary
Fagron	Xylose	16A25-B05-328261	40	not necessary
Fagron	Xylose	13A31-B40-292395	60	not necessary
Flora Apothek...	Amygdalin	13I21802	40	1405619
Flora Apothek...	Calcium colamine...	14384	40	from supplier
Flora Apothek...	Glycyrrhizic acid	10588	40	from supplier
Galke	Liquorice root, ...	22728	40	from supplier
Gatt-Koller	Metoclopramide h...	20170236	40	not necessary
Gatt-Koller	Sulfogaiacol hyd...	1156/04152115	60	not necessary
Gatt-Koller	Sulfogaiacol hyd...	4543/01170317	40	AR-17-FG-012520-01
Hedinger	Citric acid mono...	202029	40	1908675
Hedinger	Citric acid mono...	202029	40	not necessary
Inresa	Amifampridine	33-15-IV-56P2749B	40	not necessary
Inresa	Amifampridine	0631350001	60	8899
Inresa	Amitriptyline hy...	005669/17112016	40	1705099
Inresa	Amitriptyline hy...	005669/06022017	40	not necessary
Inresa	Metoprolol tartr...	309905RII	60	AR-15-FG-006635-01
Inresa	Metoprolol tartr...	405005RII	60	not necessary
Inresa	Omeprazole	OPFP14017	60	1507407
Inresa	Omeprazole	OPFP14017	60	not necessary
Inresa	Omeprazole	OPFP14017	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Inresa	Omeprazole	0PFP14017	40	not necessary
Katwijk Chemi...	Ethosuximide	451717	40	not necessary
Katwijk Chemi...	Ethosuximide	451645	40	1807457
Klenk	Hay flower	1400A130304	40	not necessary
Klenk	Potato starch	0094A170123	40	not necessary
Klenk	Tormentil	6911B141114	60	AR-15-FG-014020-01
Klenk	Wax, yellow	0290A170314	40	1908177
Klenk	Wheat starch	0104A141114	60	not necessary
Klenk	Wheat starch	0104A160512	40	not necessary
Safic-Alcan	Amygdalin, bitte...	BAE150902	60	from supplier
Safic-Alcan	Amygdalin, bitte...	BAE150902	60	from supplier
Safic-Alcan	Dimercaptosuccin...	20170210	40	from supplier
Sigma-Aldrich	Betamethasone ac...	LRAA8500	40	not necessary
Sigma-Aldrich	Betamethasone so...	MKBS9676V	40	not necessary
Sigma-Aldrich	Diclofenac potas...	D1039990	40	not necessary
Sigma-Aldrich	Fenticonazole ni...	1.3	40	not necessary
Th. Geyer	Calcium colamine...	33299	60	from supplier
Th. Geyer	Calcium colamine...	66606	40	from supplier
Th. Geyer	Mg-Ca-K colamine...	80919	40	from supplier
Th. Geyer	Mg-Ca-K colamine...	36985	60	from supplier
Th. Geyer	Mg-Ca-K colamine...	46816	40	from supplier
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	CBD-API16-004	40	2017A-03911
THC-Pharma GmbH	Cannabidiol	CBD-API16-003	40	2017B-03911
THC-Pharma GmbH	Cannabidiol	CBD-API17-003	40	not necessary
THC-Pharma GmbH	Cannabidiol	CBD-API17-003	40	not necessary
Yino Pharma	Heparin sodium	YN01011170904	40	not necessary
Yino Pharma	Heparin sodium	YN01011170502	40	not necessary
Yino Pharma	Heparin sodium	YN01011170502	40	not necessary
Yino Pharma	Heparin sodium	YN01011170904	40	not necessary
Yino Pharma	Heparin sodium	YN01011170904	40	not necessary
Yino Pharma	Sodium chondroit...	YN010511806006	40	not necessary

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Appendix B: Additional validation samples (*Type B*)

It is necessary that spectra also enter the validation which cannot be identified with this model. In this manner, it is verified that the model also rejects unknown substances. The spectra for these samples were recorded by *HiperScan GmbH*. They are allocated to the *Type B*. They also include the calibration spectra for other models.

The samples originate from 5087 batches. From these 239 317 spectra were recorded. The spectra recorded on independent samples of substances which can be identified with the model are listed respectively in the section *Type B* for the individual substances and do not appear again elsewhere in this list.

Supplier	Substance	Batch	Spectra	Certificate
abcr GmbH	Amygdalin	1310927	40	from supplier
abcr GmbH	Amygdalin	1310579	40	from supplier
abcr GmbH	Dimercaptosuccin...	1005874	39	from supplier
abcr GmbH	Glycyrrhizic acid	1222024	40	from supplier
AppliChem	Papain	1G000215	40	1407433
Audor	Alfatradiol	AP33-16-I-62P2985A	40	not necessary
Audor	Alfatradiol	AP33-18-I-68P3620A	40	not necessary
Audor	Chlorprothixene ...	AP05061117C	40	1810394
Audor	Chlorprothixene ...	AP04120418C	40	not necessary
Audor	Cyproterone acet...	APCPA140201CP	20	not necessary
Audor	Cyproterone acet...	APX03-181002-m0C	40	not necessary
Audor	Dithranol	APOP_DNL/11/17/001D	40	not necessary
Audor	Fluconazole	AP061502911F	40	not necessary
Audor	Fluconazole	AP061502911F	60	not necessary
Audor	Fluconazole	AP061502911F	60	AR-16-FG-008733-01
Audor	Fluconazole	APAFINC0217023F	40	not necessary
Audor	Fluconazole	AP061502911F	40	AR-16-FG-014225-01
Audor	Fluconazole	AP061502911F	40	not necessary
Audor	Fluconazole	AP061502911F	40	not necessary
Audor	Prasterone	APZDHA16003-1MD	40	not necessary
Audor	Tropicamide	APFP/TD003/09/17/TR	40	not necessary
Audor	Uridine	APIB1600169AU	40	not necessary
Audor	Uridine	APIB160018211AU	40	not necessary
Audor	Uridine	APIB1600169AU	40	not necessary
Auropure Life...	Amygdalin, bitte...	PA02099-160822	30	from supplier
Bombastus	Bitter fennel	289870	60	AR-15-FG-008161-01
Bombastus	Hard fat in past...	306363	6	not necessary
Bombastus	Hard fat in past...	18006248	40	not necessary
Bombastus	Maize style, cut	274920	40	not necessary
Caelo	Acetylcysteine	163340	40	20170208*
Caelo	Acetylsalicylic ...	173098	40	20171120*
Caelo	Acorus root, nat...	12076404	40	not necessary
Caelo	Acorus root, pee...	12295902	40	not necessary
Caelo	Acorus root, pee...	12239001	40	not necessary
Caelo	Agarikon	14190702	60	not necessary
Caelo	Agarikon	11024502	40	not necessary
Caelo	Agarikon	13130602	48	1405003
Caelo	Agarikon	14190706	60	not necessary
Caelo	Agnus castus fruit	13193401	40	AR-15-FG-006453-01
Caelo	Agnus castus fruit	12274003	40	1501649
Caelo	Agnus castus fruit	15318705	60	pending
Caelo	Agrimony, cut	11374004	40	not necessary
Caelo	Agrimony, cut	10325505	45	not necessary
Caelo	Alchemilla alpin...	12147201	40	not necessary
Caelo	Alchemilla, cut	12153306	40	not necessary
Caelo	Alfatradiol	161804	40	20160624*
Caelo	Alfatradiol	191598	40	20190812*
Caelo	Alfatradiol	191599	40	20190812*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Alfatradiol	191598	40	20190812*
Caelo	Alfatradiol	13201601	60	not necessary
Caelo	Alfatradiol	161804	40	20160624*
Caelo	Allantoin	181586	40	20180625*
Caelo	Allspice, powder...	12042703	40	not necessary
Caelo	Almonds, bitter,...	13220401	40	1406296
Caelo	Aloes, cape, cut	14087106	60	pending
Caelo	Aloes, cape, cut	11048808	40	pending
Caelo	Aloes, cape, pow...	11132115	40	1404400
Caelo	Aloes, cape, pow...	11132115	40	not necessary
Caelo	Aloes, cape, pow...	12183008	40	not necessary
Caelo	Aluminium oxide,...	12163019	40	not necessary
Caelo	Aluminium oxide,...	13119121	40	not necessary
Caelo	Aluminium potass...	12144605	40	from supplier
Caelo	Ambroxol hydroch...	190618	40	20190321*
Caelo	Ambroxol hydroch...	172860	40	20171103*
Caelo	Ammonium chloride	12356101	40	not necessary
Caelo	Ammonium hydroge...	182108	40	20180926*
Caelo	Angelica root	12080404	40	1403053
Caelo	Angelica root	12168701	40	pending
Caelo	Angelica root	12080409	40	not necessary
Caelo	Aniseed	11370301	40	not necessary
Caelo	Aniseed	13333903	40	not necessary
Caelo	Aniseed, pounded	12111301	40	not necessary
Caelo	Aniseed, pounded	13357101	40	not necessary
Caelo	Aniseed, powdered	11306203	40	not necessary
Caelo	Apple, cut	11229807	60	from supplier
Caelo	Arnica flower	11285702	40	not necessary
Caelo	Arnica flower	11259801	40	not necessary
Caelo	Artichoke leaf, ...	12173302	40	not necessary
Caelo	Artichoke leaf, ...	13310707	40	not necessary
Caelo	Ascorbic acid	190262	40	20190424*
Caelo	Ash leaf, cut	11031303	40	not necessary
Caelo	Asparagus root, ...	10299010	40	pending
Caelo	Atropine sulfate	170264	40	20170324*
Caelo	Atropine sulfate	172732	30	20171020*
Caelo	Atropine sulfate	16326802	40	not necessary
Caelo	Atropine sulfate	172732	40	20171020*
Caelo	Atropine sulfate	170264	40	20170324*
Caelo	Atropine sulfate	16326802	40	not necessary
Caelo	Barberries	12217109	60	pending
Caelo	Barberries	11058309	40	pending
Caelo	Barberries	14057312	60	1508391
Caelo	Barberry root ba...	12112601	40	not necessary
Caelo	Barium sulfate	171113	40	20170509*
Caelo	Basil, rubbed	12066404	50	not necessary
Caelo	Bean shells, cut	11030601	50	not necessary
Caelo	Bearberry leaf, ...	12186502	40	not necessary
Caelo	Bentonite	12237603	40	not necessary
Caelo	Benzalkonium chl...	16171603	40	not necessary
Caelo	Benzocaine	190220	40	20190322*
Caelo	Benzocaine	182778	40	20181204*
Caelo	Benzocaine	172174	40	20170928*
Caelo	Benzocaine	161740	40	20160809*
Caelo	Benzoic acid	162701	40	20161013*
Caelo	Benzoic acid	190752	40	20190502*
Caelo	Betamethasone di...	173536	25	20180111*
Caelo	Betamethasone di...	173536	40	20180111*
Caelo	Betamethasone di...	13236803	60	not necessary
Caelo	Betamethasone di...	171059	20	20170516*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Betamethasone di...	171059	40	20170516*
Caelo	Betamethasone va...	172716	40	20171017*
Caelo	Betamethasone va...	172716	20	20171017*
Caelo	Betamethasone va...	13225507	6	1706018
Caelo	Betamethasone va...	170787	20	20170418*
Caelo	Betamethasone va...	170787	40	20170418*
Caelo	Betony herb, cut	12007205	60	not necessary
Caelo	Betony herb, cut	12007205	40	not necessary
Caelo	Bifonazole	182771	40	20181114*
Caelo	Bilberry leaf, cut	12140206	40	not necessary
Caelo	Biotin	183009	40	20190107*
Caelo	Biotin	173152	40	20171129*
Caelo	Birch leaf, cut	12193105	40	not necessary
Caelo	Bismuth gallath,...	191063	40	20190524*
Caelo	Bitter fennel	161603	60	20160718*
Caelo	Bitter fennel	12167102	40	not necessary
Caelo	Bitter fennel	12094708	40	pending
Caelo	Bitter fennel	17020610	40	1708048
Caelo	Bitter fennel	13413702	40	1412549
Caelo	Bitter fennel	14376403	60	not necessary
Caelo	Bitter orange fl...	12266605	40	not necessary
Caelo	Bitter orange fl...	12266703	40	not necessary
Caelo	Bitter orange pe...	12067405	40	not necessary
Caelo	Black cohosh, cut	13237903	40	not necessary
Caelo	Black cohosh, cut	10281112	40	not necessary
Caelo	Black cohosh, cut	13237901	40	pending
Caelo	Black cumin seed	11233503	40	1403495
Caelo	Black cumin seed	161157	60	not necessary
Caelo	Black mustard seed	12189303	40	1403472
Caelo	Black mustard seed	15120405	60	AR-17-FG-007964-01
Caelo	Black mustard se...	13358304	40	pending
Caelo	Black mustard se...	16051807	60	9887
Caelo	Black pepper	11193905	40	not necessary
Caelo	Black pepper	12192311	40	1403490
Caelo	Black pepper	13221114	40	not necessary
Caelo	Blackberry leaf,...	11030707	40	not necessary
Caelo	Blackcurrant lea...	12168804	40	not necessary
Caelo	Blackcurrants	12272801	40	1403468
Caelo	Blackcurrants	14212509	60	20140829
Caelo	Blackcurrants	14212502	60	not necessary
Caelo	Blackthorn flower	12178101	40	not necessary
Caelo	Blessed thistle,...	12038203	40	not necessary
Caelo	Blue bottle flow...	13241902	60	AR-15-FG-005838-01
Caelo	Blue bottle flow...	10291105	40	pending
Caelo	Bogbean leaf, cut	12226101	40	not necessary
Caelo	Boldo leaf, cut	13367904	40	not necessary
Caelo	Borage herb, cut	11206712	40	not necessary
Caelo	Boric acid	15242903	40	not necessary
Caelo	Boric acid	191943	40	20190806*
Caelo	Bucho leaf	12034903	40	not necessary
Caelo	Budesonide, micr...	180749	40	20180306*
Caelo	Budesonide, micr...	191381	40	20190808*
Caelo	Budesonide, micr...	170788	40	20170512*
Caelo	Budesonide, micr...	191381	24	20190808*
Caelo	Budesonide, micr...	170788	20	20170512*
Caelo	Budesonide, micr...	15344004	7	not necessary
Caelo	Budesonide, micr...	11351706	40	not necessary
Caelo	Budesonide, micr...	180749	20	20180306*
Caelo	Buprenorphine hy...	16143001	40	1919M-07774
Caelo	Buprenorphine hy...	16143001	40	2017A-07774

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Buprenorphine hy...	16143001	40	1919M-07774
Caelo	Buprenorphine hy...	16143001	40	0718E-07774
Caelo	Buprenorphine hy...	16143001	40	2017A-07774
Caelo	Buprenorphine hy...	16143001	40	0718E-07774
Caelo	Buprenorphine hy...	12172908	60	0814I-07774
Caelo	Buprenorphine hy...	12172908	40	0814I-07774
Caelo	Burdock root, cut	11100801	40	not necessary
Caelo	Burnet root	12027410	40	not necessary
Caelo	Burnet root	12024406	40	not necessary
Caelo	Burnet root	12027407	40	not necessary
Caelo	Butylated hydrox...	180291	40	20180301*
Caelo	Butylated hydrox...	163163	40	20161123*
Caelo	Caffeine	15082303	40	not necessary
Caelo	Calcium ascorbate	183226	40	20190111*
Caelo	Calcium chloride...	182411	35	20181022*
Caelo	Calcium citrate ...	181608	40	20180608*
Caelo	Calcium gluconate	154356	40	20160225*
Caelo	Calcium gluconate	181010	40	20180514*
Caelo	Calcium gluconate	190240	40	20190221*
Caelo	Calcium glycerop...	12012417	40	not necessary
Caelo	Calcium glycerop...	12012411	40	not necessary
Caelo	Calcium glycerop...	12012413	40	not necessary
Caelo	Calcium glycerop...	12012417	40	not necessary
Caelo	Calcium glycerop...	162310	40	20160819*
Caelo	Calcium hydrogen...	181117	40	20180528*
Caelo	Calcium lactate ...	180734	40	20180309*
Caelo	Calcium lactate ...	171377	40	20170630*
Caelo	Calcium oxide	182244	39	20181106*
Caelo	Calcium pantothe...	15344217	40	not necessary
Caelo	Calcium pantothe...	162516	40	20160901*
Caelo	Calcium sulfate ...	163363	40	20161214*
Caelo	Calcium sulfate ...	153451	30	20151103*
Caelo	Calendula flower	11221404	40	not necessary
Caelo	Calendula flower...	11187605	40	not necessary
Caelo	Campeche wood	11074105	40	from supplier
Caelo	Camphor, racemic	180198	40	20180131*
Caelo	Camphor, racemic	172883	40	20171206*
Caelo	Camphor, racemic	183279	40	20190104*
Caelo	Camphor, racemic	16209309	40	not necessary
Caelo	Captopril	170772	40	20170327*
Caelo	Captopril	170772	40	20170327*
Caelo	Captopril	190118	40	20190312*
Caelo	Captopril	15189502	60	AR-15-FG-014855-01
Caelo	Captopril	172102	40	20170809*
Caelo	Caraway fruit	13306006	60	AR-16-FG-000540-03
Caelo	Caraway fruit	13397104	40	AR-15-FG-001526-01
Caelo	Caraway fruit	11003217	40	not necessary
Caelo	Caraway fruit	15097115	60	AR-17-FG-006649-01
Caelo	Caraway fruit	13397104	60	AR-15-FG-001527-01
Caelo	Caraway fruit, p...	12116904	40	not necessary
Caelo	Caraway fruit, p...	16073201	60	AR-17-FG-006648-01
Caelo	Caraway fruit, p...	13367101	40	not necessary
Caelo	Caraway fruit, p...	14358801	60	AR-15-FG-014914-01
Caelo	Carbachol	11150302	31	1401147
Caelo	Carbomers 35000 ...	173614	40	20180110*
Caelo	Carbomers 35000 ...	190245	40	20190214*
Caelo	Carbomers 35000 ...	180924	40	20180404*
Caelo	Carbomers 35000 ...	172520	40	20171020*
Caelo	Carbomers 35000 ...	183074	40	20190130*
Caelo	Cardamom fruit	11309107	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Cardamom fruit	13260504	40	not necessary
Caelo	Carmellose sodiu...	172446	40	20171016*
Caelo	Carmellose sodiu...	17062304	40	not necessary
Caelo	Carnauba wax	183330	40	20190114*
Caelo	Carnauba wax	170903	40	20170512*
Caelo	Castoreum, powde...	13410403	40	not necessary
Caelo	Castoreum, powde...	13410403	40	not necessary
Caelo	Castoreum, powde...	12137409	40	1402590
Caelo	Celery root, cut	10270002	40	not necessary
Caelo	Centaurry, cut	11097602	40	not necessary
Caelo	Cetomacrogol emu...	160804	40	20160323*
Caelo	Cetomacrogol emu...	160804	40	20160323*
Caelo	Cetomacrogol emu...	160804	40	20160323*
Caelo	Cetostearyl alco...	170492	40	20170302*
Caelo	Cetostearyl alco...	172979	40	20171121*
Caelo	Cetostearyl alco...	190552	40	20190404*
Caelo	Cetyl alcohol	170473	40	20170302*
Caelo	Cetyl alcohol	173155	40	20171121*
Caelo	Cetyl palmitate 15	13295919	40	not necessary
Caelo	Cetyl palmitate 15	181499	40	20180621*
Caelo	Cetyl palmitate 15	182737	40	20181221*
Caelo	Cetyl palmitate 15	13295919	20	pending
Caelo	Chickweed, cut	13254411	40	1412550
Caelo	Chickweed, cut	13040502	40	1403482
Caelo	Chicory herb, cut	12034303	60	not necessary
Caelo	Chicory herb, cut	12354302	100	not necessary
Caelo	Chloral hydrate	172506	40	20171004*
Caelo	Chloral hydrate	172477	40	20170928*
Caelo	Chloral hydrate	154483	40	20160106*
Caelo	Chloral hydrate	190914	40	20190425*
Caelo	Chloral hydrate	172477	40	20170928*
Caelo	Chloramphenicol	180178	40	20180308*
Caelo	Chloramphenicol	162786	40	20161006*
Caelo	Chloramphenicol	15358906	40	not necessary
Caelo	Chlorhexidine di...	15438914	40	not necessary
Caelo	Chlorhexidine di...	152610	40	20150805*
Caelo	Chlorhexidine di...	171858	40	20170628*
Caelo	Chlorocresol	181360	40	20180615*
Caelo	Chlorocresol	15114107	40	not necessary
Caelo	Chlorocresol	1211107	40	not necessary
Caelo	Cholesterol	13148703	40	not necessary
Caelo	Cholesterol	190549	40	20190416*
Caelo	Choline chloride	181082	40	20180716*
Caelo	Choline chloride	172053	40	20170829*
Caelo	Ciclopirox olamine	12040906	40	not necessary
Caelo	Ciclopirox olamine	182102	40	20180928*
Caelo	Ciclosporin	171583	40	20170623*
Caelo	Ciclosporin	173327	5	20180111*
Caelo	Ciclosporin	15278907	40	not necessary
Caelo	Ciclosporin	12317707	30	not necessary
Caelo	Cinchona bark, cut	11238512	40	1403463
Caelo	Cinchona bark, cut	14009304	40	AR-15-FG-007677-01
Caelo	Cinchona bark, cut	15301903	60	AR-17-FG-007831-01
Caelo	Cinchona bark, cut	13031007	40	not necessary
Caelo	Cinnamon rods (C...	13241807	40	not necessary
Caelo	Cinnamon, cut	16114401	60	1704441
Caelo	Cinnamon, cut	13304414	40	not necessary
Caelo	Cinnamon, cut	14175010	60	1505349
Caelo	Cinnamon, cut	12309902	40	pending
Caelo	Citric acid	180976	40	20180327*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Citric acid	183121	40	20181220*
Caelo	Citric acid	171626	40	20170628*
Caelo	Citric acid	181292	40	20180606*
Caelo	Citric acid	183071	40	20181220*
Caelo	Citric acid	171627	40	20170614*
Caelo	Citric acid	162086	40	20160720*
Caelo	Citric acid	162154	40	20160725*
Caelo	Citric acid	192036	40	20190808*
Caelo	Citric acid	171627	40	20170614*
Caelo	Citric acid	160387	60	20160316*
Caelo	Citric acid	192073	40	20190808*
Caelo	Citric acid	171626	40	20170628*
Caelo	Citric acid mono...	161865	40	20160706*
Caelo	Citric acid mono...	161865	40	20160706*
Caelo	Clindamycin hydr...	172093	40	20170818*
Caelo	Clindamycin hydr...	153928	40	20151204*
Caelo	Clindamycin-2-di...	172964	40	20171109*
Caelo	Clindamycin-2-di...	163265	40	20161214*
Caelo	Clindamycin-2-di...	15435801	60	not necessary
Caelo	Clindamycin-2-di...	14257305	60	AR-15-FG-015412-01
Caelo	Clindamycin-2-di...	14257305	60	AR-15-FG-015412-01
Caelo	Clioquinol	172046	40	20170904*
Caelo	Clioquinol	16307206	40	not necessary
Caelo	Clioquinol	182782	40	20190103*
Caelo	Clioquinol	13114413	40	not necessary
Caelo	Clobetasol propi...	170967	40	20170524*
Caelo	Clobetasol propi...	161067	60	20160422*
Caelo	Clobetasol propi...	161067	48	20160422*
Caelo	Clobetasol propi...	170967	20	20170524*
Caelo	Clotrimazole	171116	40	20170522*
Caelo	Clotrimazole	190162	40	20190220*
Caelo	Clotrimazole	17111616	40	not necessary
Caelo	Clotrimazole	12320815	30	not necessary
Caelo	Clotrimazole	190161	40	20190220*
Caelo	Clove	12095805	40	pending
Caelo	Clove	16207804	60	1708676
Caelo	Cocoa butter	173149	40	20171211*
Caelo	Cocoa butter	181351	40	20180514*
Caelo	Cocoa butter	182544	40	20181016*
Caelo	Cocoa butter	17055907	7	not necessary
Caelo	Colophony	173014	40	20180102*
Caelo	Colophony	173014	40	20180102*
Caelo	Common wormwood ...	12016705	40	not necessary
Caelo	Copper (II) sulf...	162409	40	20160825*
Caelo	Coriander	16010212	60	AR-17-FG-015965-01
Caelo	Coriander	13299403	40	1412548
Caelo	Coriander	12080904	40	pending
Caelo	Corn starch	170474	40	20170309*
Caelo	Corn starch	16138505	40	not necessary
Caelo	Corn starch	172873	40	20171129*
Caelo	Cottonwood bud	12027505	40	not necessary
Caelo	Couch grass rhiz...	12186702	40	not necessary
Caelo	Cowberry leaf, cut	11334302	40	not necessary
Caelo	Cumin	11378308	30	not necessary
Caelo	Cyanocobalamin	13342807	30	AR-15-FG-014917-01
Caelo	Cyanocobalamin	162215	40	20160815*
Caelo	Cyanocobalamin	16015411	40	not necessary
Caelo	Cyanocobalamin	09351704	40	not necessary
Caelo	Daisy flower	11068306	40	not necessary
Caelo	Damiana leaf, cut	11026608	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Dandelion herb w...	11217109	40	not necessary
Caelo	Dandelion root, ...	13329701	40	not necessary
Caelo	Dandelion root, ...	12363102	40	1403708
Caelo	Dandelion root, ...	13329708	40	1411089
Caelo	Dapsone	16073502	40	not necessary
Caelo	Dapsone	182792	40	20181214*
Caelo	Dapsone	16073502	40	not necessary
Caelo	D-Camphor	190097	40	20190124*
Caelo	D-Camphor	162354	40	20160809*
Caelo	Dequalinium chlo...	180794	40	20180314*
Caelo	Dequalinium chlo...	12077503	40	not necessary
Caelo	Dequalinium chlo...	15218406	40	not necessary
Caelo	Dequalinium chlo...	12361906	60	not necessary
Caelo	Dequalinium chlo...	154028	40	20160219*
Caelo	Dequalinium chlo...	12361906	60	not necessary
Caelo	Dequalinium chlo...	12361906	60	not necessary
Caelo	Dequalinium chlo...	12361906	60	not necessary
Caelo	Dequalinium chlo...	12361906	60	not necessary
Caelo	Devil's claw root	13298202	40	not necessary
Caelo	Devil's claw root	12031806	40	not necessary
Caelo	Dexamethasone	172856	40	20171102*
Caelo	Dexamethasone	13084306	60	pending
Caelo	Dexamethasone	171042	24	20170424*
Caelo	Dexamethasone	172856	40	20171102*
Caelo	Dexamethasone	15054627	20	not necessary
Caelo	Dexamethasone	190865	40	20190514*
Caelo	Dexamethasone	171042	40	20170424*
Caelo	Dexamethasone	190865	40	20190514*
Caelo	Dexamethasone ac...	191192	40	20190705*
Caelo	Dexamethasone ac...	161066	60	20160426*
Caelo	Dexamethasone ac...	161066	42	20160426*
Caelo	Dexamethasone ac...	181730	40	20180712*
Caelo	Dexamethasone ac...	181730	20	20180712*
Caelo	Dexamethasone ac...	191192	40	20190705*
Caelo	Dexamethasone ac...	14200207	15	not necessary
Caelo	Dill fruit	14018705	60	pending
Caelo	Dill fruit	12289901	40	1403464
Caelo	Dill fruit	12342906	40	not necessary
Caelo	Dill herb, cut	11295107	40	pending
Caelo	Dill herb, cut	11295103	40	not necessary
Caelo	Dill herb, cut	13007602	40	1402204
Caelo	Dill herb, cut	13007610	40	not necessary
Caelo	Diltiazem hydroc...	173370	40	20180115*
Caelo	Diltiazem hydroc...	190841	40	20190507*
Caelo	Diltiazem hydroc...	170432	40	20170331*
Caelo	Diltiazem hydroc...	183011	40	20181218*
Caelo	Diphenhydramine ...	183010	40	20181213*
Caelo	Diphenhydramine ...	172276	40	20170919*
Caelo	Diphenhydramine ...	15072617	40	not necessary
Caelo	Doxylamine hydro...	181277	40	20180613*
Caelo	Doxylamine hydro...	172047	40	20170830*
Caelo	Dry beads Caelo...	11005804	40	from supplier
Caelo	Dry beads Caelo...	12382203	40	from supplier
Caelo	Echinacea pallid...	12244501	40	not necessary
Caelo	Echinacea pallid...	14120602	40	not necessary
Caelo	Elder berries	12267807	40	1403466
Caelo	Elder berries	13367704	40	1410179
Caelo	Elder flower, ru...	12221902	40	not necessary
Caelo	Elecampane rhizo...	12113303	40	not necessary
Caelo	Eleutherococcus,...	10381814	30	not necessary
Caelo	Elm bark, cut	10347205	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Emser salt, arti...	181633	40	20180621*
Caelo	Emulsifying ceto...	171980	40	20170721*
Caelo	Emulsifying ceto...	180180	40	20180222*
Caelo	Ephedrine hydroc...	171736	40	20170627*
Caelo	Epinephrine hydr...	181795	40	20180816*
Caelo	Epinephrine hydr...	171115	40	20170510*
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	12337805	40	not necessary
Caelo	Epinephrine hydr...	14007802	60	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	162764	40	20161103*
Caelo	Epinephrine hydr...	171115	40	20170510*
Caelo	Epinephrine hydr...	171115	40	20170510*
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	171115	40	20170510*
Caelo	Epinephrine hydr...	162764	40	20161103*
Caelo	Epinephrine hydr...	181795	40	20180816*
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Epinephrine hydr...	14007802	20	not necessary
Caelo	Equisetum stem, ...	10375515	40	not necessary
Caelo	Erythromycin	160795	60	20160408*
Caelo	Erythromycin	160797	60	20160408*
Caelo	Erythromycin	160796	60	20160408*
Caelo	Erythromycin	13175601	48	not necessary
Caelo	Estradiol benzoate	15078608	40	not necessary
Caelo	Estradiol benzoate	13277607	40	not necessary
Caelo	Estradiol benzoate	181444	40	20180614*
Caelo	Estradiol hemihy...	13285802	60	1510622
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	60	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	170582	30	20170425*
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	181443	40	20180529*
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	170582	40	20170425*
Caelo	Estradiol hemihy...	13285802	20	not necessary
Caelo	Estradiol hemihy...	181443	40	20180529*
Caelo	Estradiol valerate	180674	40	20180316*
Caelo	Estradiol valerate	13264001	16	not necessary
Caelo	Estradiol valerate	180674	40	20180316*
Caelo	Estradiol valerate	13264001	15	not necessary
Caelo	Estradiol valerate	13264001	15	not necessary
Caelo	Estradiol valerate	13264001	17	not necessary
Caelo	Estradiol valerate	13264001	14	not necessary
Caelo	Estradiol valerate	13264001	15	not necessary
Caelo	Estradiol valerate	13264001	15	not necessary
Caelo	Estriol	161382	40	20160523*
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	20	not necessary
Caelo	Estriol	13256001	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Estriol	161382	40	20160523*
Caelo	Ethacridine lact...	172476	40	20171024*
Caelo	Ethacridine lact...	172476	40	20171024*
Caelo	Ethacridine lact...	183146	40	20190102*
Caelo	Eucalyptus leaf,...	11307706	40	not necessary
Caelo	Fampridine	172170	40	20170828*
Caelo	Fampridine	190675	40	20190506*
Caelo	Fenugreek	11313204	40	not necessary
Caelo	Fenugreek, powde...	12149605	40	not necessary
Caelo	Fenugreek, powde...	15308704	60	AR-17-FG-006625-01
Caelo	Flax weed herb, ...	14005309	60	not necessary
Caelo	Flax weed herb, ...	13007505	40	1403462
Caelo	Fluocinolone aceto...	172130	12	20170809*
Caelo	Fluocinolone aceto...	172130	40	20170809*
Caelo	Frangula bark, cut	12185901	40	not necessary
Caelo	Frankincense	13112801	40	not necessary
Caelo	Frankincense	12120005	40	1402495
Caelo	Fraxinella root,...	11314705	40	pending
Caelo	Fraxinella root,...	14110802	60	not necessary
Caelo	Fraxinella root,...	10189103	40	not necessary
Caelo	Fraxinella root,...	11314705	80	not necessary
Caelo	Fructose	171439	40	20170627*
Caelo	Fructose	191699	40	20190806*
Caelo	Fructose	171439	40	20170627*
Caelo	Fructose	182916	40	20181214*
Caelo	Fructose	15271213	40	not necessary
Caelo	Fumitory, cut	12103402	40	not necessary
Caelo	Galangal root	12330308	40	not necessary
Caelo	Galangal root	12330305	40	1403491
Caelo	Galangal root	11375909	40	1403324
Caelo	Galangal root	11375904	40	not necessary
Caelo	Galangal root	12125001	40	not necessary
Caelo	Galangal root	14159504	60	1511617
Caelo	Garden violet, cut	13038703	40	1403481
Caelo	Garden violet, cut	10292002	40	pending
Caelo	Garden violet, cut	13232103	40	not necessary
Caelo	Gelatine, white	172772	40	20171215*
Caelo	Gelatine, white	154419	40	20160222*
Caelo	Gentamicin sulfate	183001	40	20190117*
Caelo	Gentamicin sulfate	170996	40	20170529*
Caelo	Gentamicin sulfate	13205003	60	not necessary
Caelo	Gentamicin sulfate	183001	40	20190117*
Caelo	Gentian root, cut	12244603	40	not necessary
Caelo	Gentian root, po...	12007106	40	not necessary
Caelo	Gentian root, po...	13410503	40	1508394
Caelo	Gentian root, po...	12007103	50	1403630
Caelo	Ginger, peeled	12180601	40	not necessary
Caelo	Ginger, peeled	13321904	40	not necessary
Caelo	Ginger, peeled	14120701	40	not necessary
Caelo	Ginkgo extract	12192505	40	1403672
Caelo	Ginkgo leaf, cut	12151703	40	not necessary
Caelo	Ginkgo leaf, cut	13340907	60	not necessary
Caelo	Ginkgo leaf, cut	13340904	40	1411088
Caelo	Ginkgo leaf, cut	16281501	60	1708047
Caelo	Ginseng root	13237805	40	not necessary
Caelo	Ginseng root	13341008	40	not necessary
Caelo	Ginseng root	11349206	40	not necessary
Caelo	Ginseng root	12080208	40	not necessary
Caelo	Glucose	16157401	40	not necessary
Caelo	Glucose	180129	40	20180118*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Glucose monohydr...	190733	40	20190411*
Caelo	Glucose monohydr...	16184807	40	not necessary
Caelo	Glycerol monoste...	191186	40	20190605*
Caelo	Glycerol monoste...	190794	40	20190418*
Caelo	Glycerol monoste...	180447	40	20180427*
Caelo	Glycerol monoste...	11272902	80	not necessary
Caelo	Glycerol monoste...	162338	40	20160822*
Caelo	Glycerol monoste...	162338	40	20160822*
Caelo	Glycerol monoste...	162338	40	20160822*
Caelo	Glycerol monoste...	162338	40	20160822*
Caelo	Glycerol monoste...	162338	40	20160822*
Caelo	Glycerol monoste...	14308310	60	not necessary
Caelo	Glycerol monoste...	14308310	60	not necessary
Caelo	Glycerol monoste...	172585	40	20170311*
Caelo	Glycerol monoste...	14308310	60	not necessary
Caelo	Glycerol monoste...	11079309	40	not necessary
Caelo	Glycine	191009	40	20190509*
Caelo	Guaiacum wood, cut	11352001	40	1406340
Caelo	Guaiacum wood, cut	13343802	40	not necessary
Caelo	Guaiacum wood, cut	11352012	40	not necessary
Caelo	Guaifenesin	190134	40	20190321*
Caelo	Hamamelis bark, ...	12361301	40	not necessary
Caelo	Hamamelis bark, ...	12361303	40	not necessary
Caelo	Hamamelis leaf, ...	13269904	40	not necessary
Caelo	Hamamelis leaf, ...	13269902	40	1405176
Caelo	Hamamelis leaf, ...	10102915	40	not necessary
Caelo	Hard fat in past...	183339	40	20190227*
Caelo	Hard fat in past...	181479	40	20180904*
Caelo	Hard fat in past...	173076	40	20180110*
Caelo	Hawthorn flower	13010403	40	not necessary
Caelo	Hawthorn flower	10157201	40	not necessary
Caelo	Hawthorn fruit	11278205	30	not necessary
Caelo	Hawthorn fruit	11183010	40	not necessary
Caelo	Hawthorn leaf wi...	13276706	40	not necessary
Caelo	Hawthorn leaf wi...	10273515	40	not necessary
Caelo	Hawthorn leaf wi...	10273513	40	not necessary
Caelo	Heather flower, ...	10270508	40	not necessary
Caelo	Heather, cut	11011209	40	not necessary
Caelo	Heavy alkaline b...	181942	40	20181010*
Caelo	Henna powder	12122302	40	from supplier
Caelo	Henna powder	12267304	40	from supplier
Caelo	Hop strobile	12008305	40	not necessary
Caelo	Hop strobile	12221802	40	not necessary
Caelo	Horse chestnut l...	11063503	40	not necessary
Caelo	Hyaluronic acid ...	172347	40	20171004*
Caelo	Hydrochlorothiaz...	14009701	40	not necessary
Caelo	Hydrocortisone	171194	20	20170530*
Caelo	Hydrocortisone	171194	30	20170530*
Caelo	Hydrocortisone	160952	40	20160415*
Caelo	Hydrocortisone	171194	40	20170530*
Caelo	Hydrocortisone	160952	40	20160415*
Caelo	Hydrocortisone	172566	40	20170928*
Caelo	Hydrocortisone	13341502	60	not necessary
Caelo	Hydrocortisone	183243	40	20190131*
Caelo	Hydrocortisone	172566	40	20170928*
Caelo	Hydrocortisone	183243	40	20190131*
Caelo	Hydrocortisone	171194	28	20170530*
Caelo	Hydrocortisone	180915	40	20180412*
Caelo	Hydrocortisone a...	182780	32	20181210*
Caelo	Hydrocortisone a...	16238706	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Hydrocortisone a...	170862	40	20170412*
Caelo	Hydrocortisone a...	16238706	20	not necessary
Caelo	Hydrocortisone a...	170862	21	20170412*
Caelo	Hydroquinone	16023105	40	not necessary
Caelo	Hydroquinone	172215	40	20170912*
Caelo	Hydroxyethylcell...	191675	40	20190716*
Caelo	Hydroxyethylcell...	181131	40	20180612*
Caelo	Hydroxyethylcell...	173196	40	20180111*
Caelo	Hydroxyethylcell...	190767	40	20190430*
Caelo	Hydroxyethylcell...	17115103	40	not necessary
Caelo	Hydroxyethylcell...	16012909	40	not necessary
Caelo	Hydroxyethylcell...	16252003	40	not necessary
Caelo	Hydroxyethylcell...	190174	40	20190411*
Caelo	Hydroxypropylcel...	15345009	40	not necessary
Caelo	Hydroxypropylcel...	17355801	40	not necessary
Caelo	Hyoscine hydrobr...	15054504	40	not necessary
Caelo	Hyoscine hydrobr...	15054504	60	not necessary
Caelo	Hyoscine hydrobr...	15054501	32	1511156
Caelo	Hyoscine hydrobr...	15054504	40	not necessary
Caelo	Hyoscine hydrobr...	15054504	30	not necessary
Caelo	Hypromellose 4000	190696	40	20190424*
Caelo	Hypromellose 4000	162244	40	20160818*
Caelo	Hypromellose 4000	171150	40	20170609*
Caelo	Ibuprofen	183003	40	20190207*
Caelo	Ibuprofen	13105108	40	not necessary
Caelo	Iceland moss, cut	13350803	40	not necessary
Caelo	Iceland moss, cut	12187106	40	1403493
Caelo	Iodine	181635	40	20180704*
Caelo	Iodoform	182495	40	20181024*
Caelo	Iodoform	180919	40	20180516*
Caelo	Iodoform	170860	30	20170707*
Caelo	Iron (II) sulfat...	11340509	40	1401568
Caelo	Iron (III) chlor...	12354703	40	1401391
Caelo	Iron (III) sulfa...	170991	40	20170523*
Caelo	Isomalt	170812	40	20170421*
Caelo	Isoxuprine hydro...	190117	40	20190221*
Caelo	Ispaghula husk	12148103	40	not necessary
Caelo	Ispaghula seed	11185807	40	not necessary
Caelo	Ivy leaf, cut	11229505	40	not necessary
Caelo	Ivy leaf, cut	13276304	40	not necessary
Caelo	Jasmine flower	12143201	40	not necessary
Caelo	Java tea, cut	13258805	40	not necessary
Caelo	Java tea, cut	12225702	40	not necessary
Caelo	Juniper	13267210	40	not necessary
Caelo	Juniper	12028609	40	not necessary
Caelo	Juniper, pounded	12278903	40	1405296
Caelo	Kaolin, heavy	11350124	40	not necessary
Caelo	Kaolin, red, pow...	12093503	40	from supplier
Caelo	Kidney vetch, co...	10249206	40	not necessary
Caelo	Knotgrass, cut	14079304	40	not necessary
Caelo	Lactose monohydr...	15236220	40	not necessary
Caelo	Lactose monohydr...	173360	40	20180110*
Caelo	Ladys bedstraw ...	11265103	40	not necessary
Caelo	Lapacho wood, cut	12239305	40	not necessary
Caelo	Lapis Pumicis	180533	40	20180222*
Caelo	Laurel leaf	12299103	40	not necessary
Caelo	Lavender flower	13336005	40	from supplier
Caelo	Lavender flower	12088003	60	not necessary
Caelo	Lemon peel, cut	12217402	40	not necessary
Caelo	Lemon verbena le...	12137505	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Lemon verbena le...	12137503	40	not necessary
Caelo	Lidocaine	171831	40	20170724*
Caelo	Lidocaine	190545	40	20190315*
Caelo	Lidocaine	162895	40	20161107*
Caelo	Lidocaine hydroc...	172048	40	20170801*
Caelo	Lidocaine hydroc...	152995	40	20151013*
Caelo	Lidocaine hydroc...	162222	40	20160811*
Caelo	Lime flower, cut	11086305	55	not necessary
Caelo	Linseed	12186902	60	pending
Caelo	Linseed	12092102	60	not necessary
Caelo	Linseed	13412801	40	not necessary
Caelo	Linseed	13296703	40	not necessary
Caelo	Liquorice root, ...	12227101	40	not necessary
Caelo	Liquorice root, ...	12227106	40	not necessary
Caelo	Liquorice	14050601	60	not necessary
Caelo	Liquorice	13068904	40	not necessary
Caelo	Liquorice	13180904	60	pending
Caelo	Liquorice root, ...	12267103	40	not necessary
Caelo	Locust bean gum,...	10376516	40	not necessary
Caelo	Locust bean gum,...	13379403	40	from supplier
Caelo	Lovage root, cut	11181609	40	not necessary
Caelo	Macrogols 1500	182290	40	20181011*
Caelo	Macrogols 1500	172221	40	20170911*
Caelo	Macrogols 4000	171190	40	20170509*
Caelo	Macrogols 4000	181974	40	20180912*
Caelo	Macrogols 4000	182154	40	20180821*
Caelo	Macrogols 4000	173039	40	20171206*
Caelo	Macrogols 4000	183153	40	20181214*
Caelo	Macrogols 4000	190093	40	20190128*
Caelo	Macrogols 4000	182448	40	20181011*
Caelo	Macrogols 4000	170686	40	20170407*
Caelo	Macrogols 4000	11172433	40	not necessary
Caelo	Macrogols 4000	191756	40	20190726*
Caelo	Macrogols 4000	172081	40	20170509*
Caelo	Magensium hydrog...	12181911	40	not necessary
Caelo	Magnesium carbon...	182282	40	20180906*
Caelo	Magnesium hydrog...	173531	40	20180112*
Caelo	Magnesium hydrog...	182746	40	20181109*
Caelo	Magnesium hydrox...	16253301	30	not necessary
Caelo	Magnesium orotat...	172157	40	20170817*
Caelo	Magnesium orotat...	172547	40	20171123*
Caelo	Magnesium oxide,...	180818	40	20180322*
Caelo	Magnesium peroxi...	190555	40	20190306*
Caelo	Magnesium stearate	183089	40	20190110*
Caelo	Magnesium stearate	162066	40	20160720*
Caelo	Magnesium stearate	181533	40	20180612*
Caelo	Magnesium sulfat...	13340814	60	not necessary
Caelo	Magnesium sulfat...	172990	40	20171129*
Caelo	Maize style, cut	13381501	40	not necessary
Caelo	Maize style, cut	10325906	40	not necessary
Caelo	Male speedwell w...	12162503	40	not necessary
Caelo	Mallow flower, cut	12187405	40	not necessary
Caelo	Mallow leaf, cut	12153403	40	not necessary
Caelo	Manganese sulfat...	171648	40	20170714*
Caelo	Manna, cut	13127007	40	pending
Caelo	Manna, cut	11308907	60	pending
Caelo	Manna, powdered	11308807	40	1403534
Caelo	Manna, powdered	14001401	40	not necessary
Caelo	Mannitol	162121	40	20160810*
Caelo	Mannitol	16014013	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Mannitol	160140	60	20160303*
Caelo	Mannitol	172588	40	20171128*
Caelo	Mannitol	162009	40	20160805*
Caelo	Mannitol	182654	40	20190110*
Caelo	Mannitol	170879	40	20170424*
Caelo	Mannitol	170450	40	20170424*
Caelo	Mannitol	160140	40	20160303*
Caelo	Mannitol	162121	30	20160810*
Caelo	Mannitol	182655	40	20181219*
Caelo	Mannitol	191654	40	20190716*
Caelo	Mannitol	170879	40	20170424*
Caelo	Marian thistle f...	11126703	40	not necessary
Caelo	Marian thistle f...	13342002	40	from supplier
Caelo	Marian thistle h...	10365601	40	not necessary
Caelo	Marshmallow leaf...	12328206	40	not necessary
Caelo	Marshmallow root...	11101905	40	not necessary
Caelo	Masterwort root,...	13341204	40	not necessary
Caelo	Masterwort root,...	12101804	60	not necessary
Caelo	Maté leaf, cut	13118004	40	not necessary
Caelo	Maté leaf, cut	11091707	40	not necessary
Caelo	Matricaria flower	12294804	40	not necessary
Caelo	Matricaria flowe...	12115805	40	not necessary
Caelo	Meadowsweet flow...	12275701	40	not necessary
Caelo	Medical coal	181885	40	20180813*
Caelo	Melatonin	190946	40	20190502*
Caelo	Melatonin	190946	40	20190502*
Caelo	Melatonin	163437	40	20170111*
Caelo	Melatonin	163437	40	20170111*
Caelo	Melissa leaf, cut	13214713	40	not necessary
Caelo	Melissa leaf, cut	12032318	40	not necessary
Caelo	Menthol	170927	40	20170418*
Caelo	Menthol	160956	60	20160412*
Caelo	Menthol	170927	40	20170418*
Caelo	Menthol	171889	40	20170719*
Caelo	Menthol	183286	40	20190219*
Caelo	Metamizol sodium...	15362815	40	not necessary
Caelo	Metamizol sodium...	170584	40	not necessary
Caelo	Metamizol sodium...	161610	40	20160714*
Caelo	Methacholine chl...	13260302	60	not necessary
Caelo	Methacholine chl...	163299	40	20161206*
Caelo	Methacholine chl...	13260302	60	pending
Caelo	Methacholine chl...	16229304	40	not necessary
Caelo	Methacholine chl...	13260302	60	not necessary
Caelo	Methacholine chl...	162293	40	20160815*
Caelo	Methacholine chl...	13260303	60	not necessary
Caelo	Methacholine chl...	163299	40	20161206*
Caelo	Methacholine chl...	163299	40	20161206*
Caelo	Methacholine chl...	16229301	40	not necessary
Caelo	Methionine	190211	40	20190306*
Caelo	Methoxsalen	191960	40	20190820*
Caelo	Methoxsalen	170781	40	20170427*
Caelo	Methoxsalen	16103014	33	not necessary
Caelo	Methoxsalen	16103014	40	not necessary
Caelo	Methoxsalen	16103003	37	not necessary
Caelo	Methoxsalen	170781	36	20170427*
Caelo	Methyl-4-hydroxy...	13268813	60	not necessary
Caelo	Methyl-4-hydroxy...	15063820	40	not necessary
Caelo	Methyl-4-hydroxy...	191302	40	20190521*
Caelo	Methyl-4-hydroxy...	171581	40	20170531*
Caelo	Methylcellulose	183318	40	20190402*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Methylcellulose	171704	40	20170720*
Caelo	Methylrosanilini...	12276308	40	1509299
Caelo	Metoclopramide h...	15421102	40	not necessary
Caelo	Metronidazole	190102	40	20190225*
Caelo	Metronidazole	13095005	40	not necessary
Caelo	Metronidazole	172905	40	20171218*
Caelo	Miconazole nitrate	161447	40	20160617*
Caelo	Miconazole nitrate	182821	40	20181114*
Caelo	Miconazole nitrate	172175	30	20170911*
Caelo	Miconazole nitrate	13185601	60	not necessary
Caelo	Minoxidil	191957	40	20190828*
Caelo	Minoxidil	172278	40	20170919*
Caelo	Mistletoe, cut	11294003	30	not necessary
Caelo	Mometasone furoa...	183341	40	20190205*
Caelo	Mometasone furoa...	163532	23	20170104*
Caelo	Mometasone furoa...	13304605	60	not necessary
Caelo	Mometasone furoa...	13304605	38	not necessary
Caelo	Mometasone furoa...	183341	21	20190205*
Caelo	Mometasone furoa...	160286	30	20160316*
Caelo	Mullein flower	11366403	40	not necessary
Caelo	Mullein flower, ...	11233402	60	not necessary
Caelo	Mupirocin calcium	161608	40	20160609*
Caelo	Myrrh	12269901	40	1403470
Caelo	Myrrh	11277111	60	1403081
Caelo	Myrrh	13108701	40	AR-15-FG-000946-01
Caelo	Naloxone hydroch...	190088	40	20190212*
Caelo	Naloxone hydroch...	15322517	40	not necessary
Caelo	Naphazoline hydr...	161741	40	20160628*
Caelo	Neomycin sulfate	13062923	45	not necessary
Caelo	Neomycin sulfate	15237602	40	not necessary
Caelo	Neostigmine brom...	17131901	40	not necessary
Caelo	Neostigmine brom...	171319	40	20170621*
Caelo	Neostigmine brom...	171319	40	20170621*
Caelo	Neostigmine brom...	17131901	40	not necessary
Caelo	Nettle leaf, cut	12184702	40	not necessary
Caelo	Nettle root, cut	11254903	40	not necessary
Caelo	Nettle seed	12131507	40	1402203
Caelo	Nettle seed	11126211	40	pending
Caelo	Nettle seed	13387110	40	1410300
Caelo	Nicotinamide	190842	40	20190405*
Caelo	Nicotinamide	161136	40	20160518*
Caelo	Nicotinamide	161136	60	20160518*
Caelo	Nicotinamide	173419	40	20180111*
Caelo	Nicotinamide	11001809	40	not necessary
Caelo	Non-ionic emulsi...	171836	40	20170628*
Caelo	Non-ionic emulsi...	162231	40	20160802*
Caelo	Nutmeg	12067304	40	not necessary
Caelo	Nutmeg	13266002	40	not necessary
Caelo	Oak bark, cut	10348617	40	not necessary
Caelo	Oak bark, cut	13360306	40	not necessary
Caelo	Oak bark, cut	11351207	40	1403492
Caelo	Oak bark, powdered	11235004	40	not necessary
Caelo	Oat straw, cut	11280008	40	not necessary
Caelo	Octenidine dihyd...	160330	60	20160316*
Caelo	Octenidine dihyd...	180174	40	20180312*
Caelo	Octenidine dihyd...	15275502	60	not necessary
Caelo	Octenidine dihyd...	163355	40	20170403*
Caelo	Octenidine dihyd...	180174	40	20180312*
Caelo	Olive leaf	12316808	60	AR-15-FG-001051-01
Caelo	Olive leaf	11270106	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Orris root	12206802	40	not necessary
Caelo	Oxalic acid	13353905	60	not necessary
Caelo	Oxetacaine	17249805	50	not necessary
Caelo	Oxytetracycline ...	173168	40	20171218*
Caelo	Pancreas powder	173226	40	20180125*
Caelo	Papaverine hydro...	163159	40	20161215*
Caelo	Papaverine hydro...	16315911	40	not necessary
Caelo	Paracetamol	170284	40	20170215*
Caelo	Paracetamol	180830	40	20180522*
Caelo	Paraffin, hard	16277704	40	not necessary
Caelo	Paraffin, microc...	15412801	40	not necessary
Caelo	Parsley root, cut	12080301	40	not necessary
Caelo	Patent blue V	16306902	40	not necessary
Caelo	Patent blue V	16306902	40	not necessary
Caelo	Peppermint leaf,...	11347607	40	not necessary
Caelo	Peppermint leaf,...	12321503	130	not necessary
Caelo	Peppermint leaf,...	13288303	40	not necessary
Caelo	Peppermint leaf,...	11294401	40	not necessary
Caelo	Phenazone	170836	40	20170511*
Caelo	Phenazone	182565	40	20181203*
Caelo	Pheniramine male...	173167	40	20171229*
Caelo	Pilocarpine hydr...	171595	40	20170606*
Caelo	Pilocarpine hydr...	171242	40	20170607*
Caelo	Pilocarpine hydr...	12217319	60	not necessary
Caelo	Pilocarpine hydr...	17124208	40	not necessary
Caelo	Pilocarpine hydr...	191220	40	20190617*
Caelo	Pilocarpine hydr...	171242	40	20170607*
Caelo	Pilocarpine hydr...	191220	40	20190617*
Caelo	Pine sprouts, cut	13256105	40	not necessary
Caelo	Pine sprouts, cut	11164708	60	not necessary
Caelo	Piper cubeba	11081903	40	not necessary
Caelo	Podophyllum resi...	12287810	40	1405242
Caelo	Poloxamers 407	181607	40	20180713*
Caelo	Poloxamers 407	182636	40	20181220*
Caelo	Poloxamers 407	172921	40	20171120*
Caelo	Poloxamers 407	170486	40	20170223*
Caelo	Polymyxin B sulf...	191131	40	20190801*
Caelo	Potassium bromide	11208211	40	not necessary
Caelo	Potassium carbon...	180995	40	20180412*
Caelo	Potassium citrate	170929	40	20170517*
Caelo	Potassium citrate	191804	40	20190729*
Caelo	Potassium citrate	180385	40	20180321*
Caelo	Potassium hydrog...	172493	40	20170911*
Caelo	Potassium hydrog...	17065603	40	not necessary
Caelo	Potassium hydrog...	11345805	40	not necessary
Caelo	Potassium hydrog...	11345805	40	not necessary
Caelo	Potassium hydrog...	17065604	40	not necessary
Caelo	Potassium hydrox...	180083	40	20180118*
Caelo	Potassium iodide	11090222	40	not necessary
Caelo	Potassium metabi...	12325906	40	not necessary
Caelo	Potassium monohy...	12366501	40	10401307
Caelo	Potassium nitrate	12236703	40	not necessary
Caelo	Potassium sodium...	172577	40	20171128*
Caelo	Potassium sodium...	15375607	40	not necessary
Caelo	Potassium sorbate	171700	40	20170711*
Caelo	Potassium sorbate	173451	40	20180109*
Caelo	Potassium sorbate	191072	40	20190503*
Caelo	Potassium sorbate	181650	40	20180716*
Caelo	Potassium tartra...	182203	40	20180907*
Caelo	Potassium tartra...	13331101	40	1403536

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Potato starch	15062921	40	not necessary
Caelo	Potato starch	173003	40	20171115*
Caelo	Povidone 25, sol...	180766	40	20180327*
Caelo	Povidone 25, sol...	182383	40	20181022*
Caelo	Povidone 25, sol...	183000	40	20190104*
Caelo	Povidone 25, sol...	172091	40	20170823*
Caelo	Povidone 25, sol...	181820	40	20180718*
Caelo	Povidone iodine	180763	40	20180503*
Caelo	Prednicarbate, m...	172204	40	20170912*
Caelo	Prednicarbate, m...	170138	40	20170202*
Caelo	Prednicarbate, m...	170138	20	20170202*
Caelo	Prednicarbate, m...	173156	20	20171129*
Caelo	Prednicarbate, m...	160874	60	20160428*
Caelo	Prednicarbate, m...	172204	20	20170912*
Caelo	Prednicarbate, m...	173156	40	20171129*
Caelo	Prednicarbate, m...	160874	30	20160428*
Caelo	Prednisolon dihy...	13304704	60	not necessary
Caelo	Prednisolone ace...	13262605	44	not necessary
Caelo	Prednisolone ace...	180673	40	20180326*
Caelo	Prednisolone ace...	180673	80	20180326*
Caelo	Prednisolone, mi...	13343105	60	not necessary
Caelo	Prednisolone, mi...	172777	38	20171212*
Caelo	Prednisolone, mi...	172777	40	20171212*
Caelo	Prednisolone, mi...	161343	3	20160601*
Caelo	Prednisone	12232716	60	not necessary
Caelo	Prednisone	171665	30	20170606*
Caelo	Prednisone	17166504	30	not necessary
Caelo	Prednisone	171665	20	20170606*
Caelo	Prednisone	12232716	60	not necessary
Caelo	Prednisone	17166504	30	not necessary
Caelo	Prilocaine hydro...	15390101	33	not necessary
Caelo	Prilocaine hydro...	15230305	60	not necessary
Caelo	Prilocaine hydro...	15390101	40	not necessary
Caelo	Primrose flower	11206806	100	not necessary
Caelo	Primrose flower	11268002	40	not necessary
Caelo	Primrose root, cut	12064506	40	not necessary
Caelo	Procaine hydroch...	161378	40	20160606*
Caelo	Procaine hydroch...	190440	40	20190305*
Caelo	Procaine hydroch...	171675	40	20170703*
Caelo	Procaine hydroch...	182764	40	20181205*
Caelo	Progesterone, mi...	15319601	40	not necessary
Caelo	Progesterone, mi...	170609	40	20170315*
Caelo	Progesterone, mi...	172778	40	20171110*
Caelo	Progesterone, mi...	172778	30	20171110*
Caelo	Progesterone, mi...	170609	20	20170315*
Caelo	Progesterone, mi...	170609	40	20170315*
Caelo	Progesterone, mi...	191287	40	20190621*
Caelo	Progesterone, mi...	17060903	20	not necessary
Caelo	Progesterone, mi...	191287	40	20190621*
Caelo	Progesterone, mi...	170609	23	20170315*
Caelo	Promethazine hyd...	170586	40	20170310*
Caelo	Promethazine hyd...	17058601	40	not necessary
Caelo	Propolis extract...	13396103	40	not necessary
Caelo	Propolis extract...	13396103	40	not necessary
Caelo	Propolis extract...	11124120	40	1403460
Caelo	Propranolol hydr...	171722	40	20170712*
Caelo	Propranolol hydr...	173330	40	20180109*
Caelo	Propyl-4-hydroxy...	171970	40	20170810*
Caelo	Propyphenazone	15348908	40	not necessary
Caelo	Propyphenazone	172776	40	20171024*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Propyphenazone	190190	40	20190315*
Caelo	Psyllium seed	13305903	40	pending
Caelo	Psyllium seed	10291206	40	pending
Caelo	Psyllium seed	10291205	80	not necessary
Caelo	Pumpkin seed, pe...	12249002	40	1404542
Caelo	Pumpkin seed, pe...	12249006	40	not necessary
Caelo	Pyridoxine hydro...	180287	40	20180221*
Caelo	Quinine dihydroc...	162696	40	20161021*
Caelo	Quinine dihydroc...	16269603	40	not necessary
Caelo	Quinine dihydroc...	16269603	40	1709386
Caelo	Quinine hydrochl...	180226	40	20180306*
Caelo	Quinine sulfate ...	170616	40	20170503*
Caelo	Quinine sulfate ...	180630	40	20180228*
Caelo	Quinine sulfate ...	170616	40	20170503*
Caelo	Quinoline yellow	16050101	40	not necessary
Caelo	Quinoline yellow	190734	40	20190425*
Caelo	Quinoline yellow	172277	30	20171020*
Caelo	Quinoline yellow	190734	40	20190425*
Caelo	Quinolinol sulfa...	16221809	40	not necessary
Caelo	Ramsons herb, cut	13072801	40	not necessary
Caelo	Ramsons herb, cut	12380801	40	1404538
Caelo	Raspberry leaf, ...	12362701	40	not necessary
Caelo	Red poppy petals...	12008104	40	not necessary
Caelo	Red-rose petal	11299204	40	not necessary
Caelo	Red-rose petal, ...	11304912	40	not necessary
Caelo	Redwood	12049604	40	not necessary
Caelo	Restharrow root,...	12330001	40	not necessary
Caelo	Rhubarb, cut	12211906	40	not necessary
Caelo	Rhubarb, cut	11057818	40	not necessary
Caelo	Rhubarb, cut	13235213	40	not necessary
Caelo	Rhubarb, powdered	11043818	60	1403280
Caelo	Rhubarb, powdered	16009901	60	9886
Caelo	Riboflavin	172722	40	20171011*
Caelo	Riboflavin	15330708	40	not necessary
Caelo	Ribwort plantain...	12362902	40	not necessary
Caelo	Rice starch	190168	40	20190322*
Caelo	Rice starch	162180	40	20160901*
Caelo	Rice starch	190169	40	20190322*
Caelo	Rifampicin	152446	40	20150831*
Caelo	Rifaximin	161053	60	20160607*
Caelo	Rose hips peel, ...	11301407	60	not necessary
Caelo	Rose hips with p...	11232204	40	not necessary
Caelo	Rose hips with p...	12064202	40	not necessary
Caelo	Rose mallow flow...	10244508	60	not necessary
Caelo	Roselle	11036509	60	not necessary
Caelo	Roselle	12007311	60	not necessary
Caelo	Rosemary leaf, cut	11193408	40	not necessary
Caelo	Rosemary leaf, cut	14159601	60	not necessary
Caelo	Saccharin sodium	182104	40	20181031*
Caelo	Saccharin sodium	15087704	40	not necessary
Caelo	Saccharin sodium	190964	40	20190424*
Caelo	Sage leaf, cut	11347503	40	not necessary
Caelo	Salicylic acid	182300	40	20180925*
Caelo	Salicylic acid	171333	40	20170517*
Caelo	Salicylic acid	160500	60	20160422*
Caelo	Salicylic acid	173015	40	20171128*
Caelo	Sarsaparilla roo...	10230402	60	not necessary
Caelo	Senna leaf, cut	12273106	40	not necessary
Caelo	Senna leaf, powd...	11193307	120	pending
Caelo	Senna leaf, powd...	12369908	40	1408079

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Senna pods, tinn...	12267605	40	not necessary
Caelo	Senna pods, tinn...	11154404	40	not necessary
Caelo	Senna pods, tinn...	13395106	40	not necessary
Caelo	Shepherd's purse...	11328208	40	not necessary
Caelo	Siam benzoin	12321314	40	not necessary
Caelo	Sildenafil citrate	190866	40	20190424*
Caelo	Silica, colloida...	181124	40	20180424*
Caelo	Silver nitrate	181276	40	20180716*
Caelo	Silver protein	190092	40	20190221*
Caelo	Silver, colloida...	173615	40	20180312*
Caelo	Silverweed, cut	11328105	40	not necessary
Caelo	Sodium acetate t...	183152	40	20181213*
Caelo	Sodium acetate t...	163007	40	20161114*
Caelo	Sodium ascorbate	162857	40	20161116*
Caelo	Sodium ascorbate	172235	40	20170911*
Caelo	Sodium ascorbate	180326	40	20180222*
Caelo	Sodium ascorbate	191393	40	20190617*
Caelo	Sodium benzoate	173530	40	20180109*
Caelo	Sodium benzoate	170558	40	20170323*
Caelo	Sodium benzoate	161371	60	20160527*
Caelo	Sodium benzoate	190541	40	20190313*
Caelo	Sodium bromide	171191	40	20170510*
Caelo	Sodium carbonate...	173613	40	20180115*
Caelo	Sodium carbonate...	16310205	40	not necessary
Caelo	Sodium citrate	182083	40	20180912*
Caelo	Sodium citrate	12045008	40	not necessary
Caelo	Sodium citrate	190890	40	20190417*
Caelo	Sodium citrate	162179	40	20160805*
Caelo	Sodium citrate	173236	40	20180115*
Caelo	Sodium cromoglic...	182657	40	20181210*
Caelo	Sodium cromoglic...	15224404	60	not necessary
Caelo	Sodium cromoglic...	173081	30	20171220*
Caelo	Sodium cyclamate	16119302	40	not necessary
Caelo	Sodium cyclamate	170411	40	20170216*
Caelo	Sodium dihydroge...	12198506	40	not necessary
Caelo	Sodium dodecyl s...	16229205	40	not necessary
Caelo	Sodium dodecyl s...	13411305	60	not necessary
Caelo	Sodium dodecyl s...	190187	40	20190218*
Caelo	Sodium dodecyl s...	181822	40	20180813*
Caelo	Sodium dodecyl s...	172363	40	20170922*
Caelo	Sodium edetate	172618	40	20171006*
Caelo	Sodium edetate	172618	40	20171006*
Caelo	Sodium iodide	161701	40	20160706*
Caelo	Sodium methyl-4...	15358709	40	not necessary
Caelo	Sodium monohydro...	12109002	40	not necessary
Caelo	Sodium nitrate	11045910	40	not necessary
Caelo	Sodium nitrite	171360	30	20170522*
Caelo	Sodium perborate...	11158903	40	not necessary
Caelo	Sodium salicylat...	16305301	40	not necessary
Caelo	Sodium thiosulfa...	181648	40	20180712*
Caelo	Softisan® 601	162364	40	20160823*
Caelo	Softisan® 601	162364	40	20160823*
Caelo	Softisan® 601	14099604	40	1504720
Caelo	Solanum dulcamar...	13254305	40	not necessary
Caelo	Solanum dulcamar...	10375602	40	not necessary
Caelo	Sorbic acid	182876	40	20181212*
Caelo	Sorbitan stearate	183125	40	20190116*
Caelo	Sorbitan stearate	171781	40	20170718*
Caelo	Sorbitan stearate	180432	40	20180221*
Caelo	Sorbitol	181636	40	20180807*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Sorbitol	190173	40	20190321*
Caelo	Sour cherry stem...	12042302	40	1503570
Caelo	Soy lecithin (in...	17118101	40	not necessary
Caelo	Soy lecithin (in...	173224	40	20171218*
Caelo	Spearmint leaf, ...	12133906	40	not necessary
Caelo	Spironolactone	182080	40	20180925*
Caelo	Spironolactone	16077305	40	not necessary
Caelo	Spironolactone	191979	40	20190731*
Caelo	Spironolactone	171165	40	20170524*
Caelo	Spironolactone	172432	40	20170905*
Caelo	Spray-dried arab...	16233105	40	not necessary
Caelo	Spray-dried arab...	16233105	40	not necessary
Caelo	Spray-dried arab...	162331	40	20160818*
Caelo	Star anise	11314104	40	not necessary
Caelo	Star anise, cut	13041401	40	1403709
Caelo	Star anise, cut	15329302	60	AR-16-FG-004770-01
Caelo	Star anise, cut	12226301	60	pending
Caelo	Stearic acid 50	180809	40	20180322*
Caelo	Stearic acid 50	163469	40	20170126*
Caelo	Stemless carline...	12163702	40	not necessary
Caelo	Strawberry leaf...	11117903	40	not necessary
Caelo	Sucrose	182574	40	20181009*
Caelo	Sucrose	182818	40	20181119*
Caelo	Sucrose	12146801	40	not necessary
Caelo	Sucrose	181546	40	20180605*
Caelo	Sulfadiazine	172222	30	20171019*
Caelo	Sulfadiazine	190683	40	20190430*
Caelo	Sulfamerazine	10218403	30	not necessary
Caelo	Sulfogaiacol hyd...	170919	40	20170511*
Caelo	Sulfogaiacol hyd...	170919	40	20170511*
Caelo	Sulfogaiacol hyd...	190921	40	20190430*
Caelo	Sweet woodruff h...	10375812	30	not necessary
Caelo	Talcum	11318312	80	1405338
Caelo	Tannin	16074504	40	not necessary
Caelo	Tartaric acid	13217803	60	not necessary
Caelo	Tartaric acid	181212	40	20180511*
Caelo	Testosterone	191113	40	20190507*
Caelo	Testosterone	171199	40	20170524*
Caelo	Testosterone	16311506	30	not necessary
Caelo	Testosterone	171199	40	20170524*
Caelo	Testosterone	171199	30	20170524*
Caelo	Testosterone	181478	40	20180706*
Caelo	Testosterone	181478	40	20180706*
Caelo	Testosterone	171199	28	20170524*
Caelo	Testosterone	191113	40	20190507*
Caelo	Testosterone pro...	191986	40	20190906*
Caelo	Testosterone pro...	191986	40	20190906*
Caelo	Testosterone pro...	15444709	30	not necessary
Caelo	Testosterone pro...	170963	40	20170427*
Caelo	Testosterone pro...	170963	29	20170427*
Caelo	Tetracaine hydro...	190119	40	20190219*
Caelo	Tetracycline hyd...	181986	40	20180822*
Caelo	Tetracycline hyd...	183083	40	20181214*
Caelo	Tetracycline hyd...	172858	40	20180112*
Caelo	Thalidomide	172100	40	20170821*
Caelo	Thalidomide	162703	40	20160929*
Caelo	Theophylline	171962	40	20170807*
Caelo	Theophylline	15203018	40	not necessary
Caelo	Theophylline-eth...	182738	40	20181203*
Caelo	Theophylline-eth...	172999	40	20171117*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Theriac, powdered	13407301	40	1410336
Caelo	Theriac, powdered	12289307	40	1405122
Caelo	Thiamine chlorid. ...	170976	40	20170523*
Caelo	Thiomersal	172173	30	20170919*
Caelo	Thyme, rubbed	11345004	60	not necessary
Caelo	Thymol	181424	40	20180613*
Caelo	Thymol	190342	40	20190319*
Caelo	Thymol	173407	40	20180110*
Caelo	Tonka beans	12145101	40	not necessary
Caelo	Tormentil	14360103	60	1606520
Caelo	Tormentil	12040708	40	1403487
Caelo	Tormentil	12040715	40	1403710
Caelo	Tormentil	13220502	40	1403049
Caelo	Tormentil	12013103	40	1403048
Caelo	Tosylchloramide ...	13361606	60	not necessary
Caelo	Tosylchloramide ...	162366	40	20160819*
Caelo	Tosylchloramide ...	172265	40	20170906*
Caelo	Tragacanth, powd. ...	13277804	60	1507149
Caelo	Tragacanth, powd. ...	16110802	60	AR-17-FG-007533-01
Caelo	Tragacanth, powd. ...	15379906	60	not necessary
Caelo	Tragacanth, powd. ...	15085602	60	not necessary
Caelo	Tretinoin	181673	40	20180628*
Caelo	Tretinoin	181674	40	20180628*
Caelo	Tretinoin	162915	40	20161109*
Caelo	Tretinoin	181673	20	20180628*
Caelo	Tretinoin	172731	40	20171019*
Caelo	Tretinoin	162171	40	20160804*
Caelo	Tretinoin	172731	40	20171019*
Caelo	Tretinoin	173057	40	20171122*
Caelo	Tretinoin	181674	40	20180628*
Caelo	Tretinoin	162171	40	20160804*
Caelo	Tretinoin	173057	40	20171122*
Caelo	Tretinoin	162915	40	20160911*
Caelo	Triamcinolone ac. ...	190619	40	20190411*
Caelo	Triamcinolone ac. ...	171206	20	20170517*
Caelo	Triamcinolone ac. ...	12301505	60	pending
Caelo	Triamcinolone ac. ...	171206	20	20170517*
Caelo	Triamcinolone ac. ...	171206	40	20170517*
Caelo	Triamcinolone ac. ...	190619	20	20190411*
Caelo	Triamcinolone ac. ...	181002	20	20180426*
Caelo	Triamcinolone ac. ...	160945	40	20160420*
Caelo	Triamcinolone ac. ...	171206	40	20170517*
Caelo	Triamcinolone ac. ...	181002	40	20180426*
Caelo	Triamcinolone ac. ...	160945	60	20160420*
Caelo	Tribromophenol b. ...	18144202	40	not necessary
Caelo	Triclosan	171498	40	20170529*
Caelo	Triclosan	182537	40	20181221*
Caelo	Triclosan	180814	40	20180322*
Caelo	Triclosan	160846	60	20160517*
Caelo	Triclosan	162808	40	20161010*
Caelo	Trometamol	173208	40	20171207*
Caelo	Turmeric root, cut	12136503	40	1405295
Caelo	Turmeric root, cut	13360202	40	1411363
Caelo	Turmeric root, cut	11160403	60	not necessary
Caelo	Turmeric, javane. ...	12185605	40	1403465
Caelo	Turmeric, javane. ...	13262904	40	not necessary
Caelo	Turmeric, javane. ...	12319206	40	not necessary
Caelo	Ubidecarenone	182298	40	20180928*
Caelo	Urea	18023003	40	not necessary
Caelo	Urea	190316	40	20190321*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Urea	180242	40	20180206*
Caelo	Urea	180230	40	20180216*
Caelo	Urea	170216	40	20170202*
Caelo	Valerian root, cut	12033809	40	not necessary
Caelo	Valerian root, p...	13113204	60	pending
Caelo	Valerian root, p...	12138005	40	1403809
Caelo	Valerian root, p...	12138008	40	1508496
Caelo	Vanillin	17117201	40	not necessary
Caelo	Vanillin	173088	40	20171221*
Caelo	Vanillin	171172	40	20170427*
Caelo	Verbena herb, cut	12187602	40	not necessary
Caelo	Walnut bark, cut	11325603	40	not necessary
Caelo	Walnut leaf, cut	11103007	40	1403091
Caelo	Wax, bleached	16163109	40	not necessary
Caelo	Wax, bleached	161631	40	20160601*
Caelo	Wax, bleached	182391	40	20181009*
Caelo	Wax, yellow	17011602	40	not necessary
Caelo	Wax, yellow	162750	40	20161007*
Caelo	Wheat starch	15421507	40	not necessary
Caelo	Wheat starch	181200	40	20180509*
Caelo	White deadnettle...	11371604	40	not necessary
Caelo	White mustard seed	13380501	40	not necessary
Caelo	White mustard seed	12361002	40	1403473
Caelo	White mustard seed	13064301	40	pending
Caelo	White mustard se...	12349601	60	1403281
Caelo	White mustard se...	12375705	40	not necessary
Caelo	White pepper	12294902	40	1403488
Caelo	White pepper	14010703	40	not necessary
Caelo	White pepper	12294905	40	not necessary
Caelo	Wild pansy, cut	10291007	60	not necessary
Caelo	Willow bark, cut	12185204	40	not necessary
Caelo	Willow herb, cut	11347311	40	not necessary
Caelo	Willowherb, smal...	13355904	60	not necessary
Caelo	Willowherb, smal...	10326106	40	not necessary
Caelo	Wool alcohols	181619	40	20180606*
Caelo	Wool alcohols	173160	40	20171120*
Caelo	Wormwood herb, cut	12355801	40	not necessary
Caelo	Xylitol	163174	40	20161207*
Caelo	Xylitol	15388101	40	not necessary
Caelo	Xylitol	181133	40	20180509*
Caelo	Xylometazoline h...	173457	40	20180109*
Caelo	Xylometazoline h...	170739	40	20170404*
Caelo	Xylometazoline h...	13344901	60	not necessary
Caelo	Xylometazoline h...	170739	40	20170404*
Caelo	Xylometazoline h...	173457	40	20180109*
Caelo	Xylometazoline h...	170739	40	20170404*
Caelo	Xylometazoline h...	16074211	40	not necessary
Caelo	Xylometazoline h...	170739	40	20170404*
Caelo	Xylose	171609	40	20170704*
Caelo	Xylose	162064	40	20160713*
Caelo	Xylose	190499	40	20190305*
Caelo	Xylose	13133402	40	not necessary
Caelo	Xylose	173188	40	20180111*
Caelo	Yarrow flower, r...	12267703	40	not necessary
Caelo	Yellow chaste we...	11293101	40	not necessary
Caelo	Zedoary root, cut	12125301	40	not necessary
Caelo	Zedoary root, po...	13165501	40	pending
Caelo	Zedoary root, po...	12344601	40	1403489
Caelo	Zedoary root, po...	14249001	60	not necessary
Caelo	Zinc acetate dih...	17020404	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Caelo	Zinc acetate dihydroxide	182046	40	20180827*
Caelo	Zinc acetate dihydroxide	170204	40	20172601*
Caelo	Zinc chloride	11345504	40	not necessary
Caelo	Zinc gluconate	161788	40	20160708*
Caelo	Zinc oxide	14038101	40	not necessary
Caelo	Zinc oxide-talcum	181487	40	20180702*
Caelo	Zinc stearate	182279	40	20181017*
Caelo	Zinc sulfate heptahydrate	12257802	40	not necessary
Carotech	Tocomax 15%	100912	40	from supplier
Dr. Behr	Quercetin	WE-13030084-02	40	not necessary
Dr. Behr	Yam root extract	WE14020005-05	40	pending
Dr. K. Hollbo...	Miconazole	029C064	40	not necessary
Euro OTC	Acetylcysteine	1801037	40	20180215*
Euro OTC	Acetylcysteine	1907008	40	20190830*
Euro OTC	Acetylcysteine	1809026	40	20181025*
Euro OTC	Acetylsalicylic acid	1805003	40	20180505*
Euro OTC	Acetylsalicylic acid	1703010	40	20170509*
Euro OTC	Acetylsalicylic acid	1805004	40	20180604*
Euro OTC	Alfatradiol	1902020	40	20190402*
Euro OTC	Alfatradiol	1902020	40	20190402*
Euro OTC	Allantoin	1811006	40	20181123*
Euro OTC	alpha-tocopheryl acetate	1705003	40	20170622*
Euro OTC	alpha-tocopheryl acetate	1810024	40	20181204*
Euro OTC	Ambroxol hydrochloride	1711021	40	20171208*
Euro OTC	Ambroxol hydrochloride	1806009	40	20180625*
Euro OTC	Ambroxol hydrochloride	1811039	40	20190116*
Euro OTC	Ambroxol hydrochloride	1609002	40	20161004*
Euro OTC	Ambroxol hydrochloride	1904001	40	20190425*
Euro OTC	Ambroxol hydrochloride	1702004	40	20170320*
Euro OTC	Ammonium chloride	1209013-01	40	not necessary
Euro OTC	Amphotericin B	1803036	40	20180418*
Euro OTC	Amphotericin B	1903002	40	20190404*
Euro OTC	Arginine	L1405020-01	60	not necessary
Euro OTC	Arginine	1902024	40	20190319*
Euro OTC	Arginine	1909005	40	20190913*
Euro OTC	Arginine	1806010	40	20180619*
Euro OTC	Arginine	1710010	40	20171019*
Euro OTC	Arginine hydrochloride	1711014	40	20171127*
Euro OTC	Arginine hydrochloride	1809012	40	20181015*
Euro OTC	Ascorbic acid	1710032	40	20171123*
Euro OTC	Ascorbic acid	1901044	40	20190211*
Euro OTC	Ascorbic acid	1707005	40	20170714*
Euro OTC	Ascorbic acid	1703025	40	20170323*
Euro OTC	Ascorbic acid	1708026	40	20170905*
Euro OTC	Ascorbic acid	1903008	40	20190402*
Euro OTC	Ascorbic acid	1703040	40	20170403*
Euro OTC	Ascorbic acid	1909018	40	20191011*
Euro OTC	Beclometasone dihydroxide	1203022-03	5	1509281
Euro OTC	Beclometasone dihydroxide	1203022-03	28	not necessary
Euro OTC	Betamethasone dihydroxide	1811007	40	20181122*
Euro OTC	Betamethasone dihydroxide	1610014	40	20161116*
Euro OTC	Betamethasone dihydroxide	1705020	27	20170727*
Euro OTC	Betamethasone dihydroxide	1801045	40	20180327*
Euro OTC	Betamethasone dihydroxide	1811007	30	20181122*
Euro OTC	Betamethasone dihydroxide	1705020	40	20170727*
Euro OTC	Betamethasone dihydroxide	1907034	40	20190805*
Euro OTC	Betamethasone dihydroxide	1801045	39	20180327*
Euro OTC	Betamethasone dihydroxide	1907034	40	20190805*
Euro OTC	Betamethasone dihydroxide	1610014	30	20161116*
Euro OTC	Betamethasone valerate	1606023	34	20160708*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Betamethasone va...	1606023	40	20160708*
Euro OTC	Betamethasone va...	1704019	40	20170529*
Euro OTC	Betamethasone va...	1704019	29	20170529*
Euro OTC	Betamethasone va...	1801024	40	20180214*
Euro OTC	Betamethasone va...	1610010	30	20161116*
Euro OTC	Betamethasone va...	1801024	40	20180214*
Euro OTC	Bifonazole	1510001	60	20151203*
Euro OTC	Bifonazole	1703041	40	20170612*
Euro OTC	Bifonazole	1806008	40	20180619*
Euro OTC	Biotin	1903019	40	20190506*
Euro OTC	Biotin	1708007	40	20171023*
Euro OTC	Biotin	1603017	40	20160427*
Euro OTC	Biotin	1802022	40	20180302*
Euro OTC	Biotin	1903020	40	20190506*
Euro OTC	Budesonide, micr...	1709022	20	20171023*
Euro OTC	Budesonide, micr...	1904004	40	20190515*
Euro OTC	Budesonide, micr...	1904004	40	20190515*
Euro OTC	Budesonide, micr...	1709022	40	20171023*
Euro OTC	Caffeine	1907006	40	20190718*
Euro OTC	Caffeine	1703006	40	20170322*
Euro OTC	Caffeine	1603009	40	20160331*
Euro OTC	Calcium ascorbate	1904010	40	20190521*
Euro OTC	Calcium ascorbate	1901003	40	20190212*
Euro OTC	Calcium ascorbate	1607033-01	40	not necessary
Euro OTC	Calcium ascorbate	1807032	40	20180829*
Euro OTC	Calcium carbonat...	1611008-01	40	not necessary
Euro OTC	Calcium citrate ...	1803031	40	20180326*
Euro OTC	Calcium citrate ...	1901002	40	20190125*
Euro OTC	Calcium citrate ...	1606001	40	20160609*
Euro OTC	Calcium citrate ...	1707025	40	20170907*
Euro OTC	Calcium gluconate	1606034	40	20160630*
Euro OTC	Calcium gluconate	1712007	40	20171220*
Euro OTC	Calcium gluconate	1905015	40	20190705*
Euro OTC	Calcium orotate ...	1009035-01	40	not necessary
Euro OTC	Calcium orotate ...	1411003	40	20160331*
Euro OTC	Calcium pantothe...	1803014	40	20180316*
Euro OTC	Calcium pantothe...	1809014	40	20181023*
Euro OTC	Calcium pantothe...	1902032	40	20190326*
Euro OTC	Capsaicin	1603005-01	40	not necessary
Euro OTC	Capsaicin	1603005-01	40	not necessary
Euro OTC	Capsaicin	1812015	40	20190128*
Euro OTC	Capsaicin	1806013-01	40	not necessary
Euro OTC	Capsaicin	118066-04	60	not necessary
Euro OTC	Capsaicin	1903032	40	20190426*
Euro OTC	Carbomers 35000 ...	1810019	40	20181120*
Euro OTC	Carbomers 35000 ...	1903007	40	20190329*
Euro OTC	Carbomers 35000 ...	1702012	40	20170427*
Euro OTC	Carbomers 35000 ...	1801044	40	20180208*
Euro OTC	Carbomers 35000 ...	1706015	40	20170802*
Euro OTC	Carbomers 35000 ...	1908029	40	20190909*
Euro OTC	Cellulose, micro...	1610030	40	20161125*
Euro OTC	Cellulose, micro...	1604029-01	40	not necessary
Euro OTC	Cellulose, micro...	1607042	40	20160902*
Euro OTC	Cellulose, micro...	1811011	40	20181114*
Euro OTC	Cellulose, micro...	1712009	40	20180112*
Euro OTC	Chloramphenicol	1704021	40	20170530*
Euro OTC	Chloramphenicol	1903028	40	20190503*
Euro OTC	Chlortetracyclin...	1803037	40	20180412*
Euro OTC	Chlortetracyclin...	1604015-01	40	not necessary
Euro OTC	Cholecalciferol ...	1901050	40	20190221*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Cholecalciferol ...	1812010	40	20190114*
Euro OTC	Cholecalciferol ...	1806003	40	20180618*
Euro OTC	Choline hydrogen...	L1208016-01	40	1402479
Euro OTC	Choline hydrogen...	1802008	40	20180215*
Euro OTC	Choline hydrogen...	1707023	40	20170907*
Euro OTC	Choline hydrogen...	1711022	40	20171129*
Euro OTC	Choline hydrogen...	L1208016-02	40	not necessary
Euro OTC	Chromium yeast 0...	L1406030-02	60	not necessary
Euro OTC	Chromium yeast 0...	L1303028-01	40	1403474
Euro OTC	Citrulline	1904013	40	20190508*
Euro OTC	Citrulline	1707003	40	20170727*
Euro OTC	Citrulline	1609013	40	20161006*
Euro OTC	Citrulline	1710039	40	20171124*
Euro OTC	Citrulline	1805028	40	20180613*
Euro OTC	Clindamycin hydr...	1811003	40	20181203*
Euro OTC	Clindamycin hydr...	1903009	40	20190402*
Euro OTC	Clindamycin hydr...	1401041-01	40	not necessary
Euro OTC	Clindamycin hydr...	1710013	40	20171107*
Euro OTC	Clioquinol	1612020	40	20170215*
Euro OTC	Clioquinol	1903010	40	20190425*
Euro OTC	Clobetasol propi...	1903014	40	20190409*
Euro OTC	Clobetasol propi...	1704020	20	20170613*
Euro OTC	Clobetasol propi...	1704020	40	20170613*
Euro OTC	Clobetasol propi...	1804002	40	20180426*
Euro OTC	Clobetasol propi...	1804002	22	20180426*
Euro OTC	Clobetasol propi...	1903014	40	20190409*
Euro OTC	Clotrimazole	1902004	40	20190318*
Euro OTC	Clotrimazole	1703028	40	20170412*
Euro OTC	Clotrimazole	1607046	40	20160823*
Euro OTC	Clotrimazole	1808028	40	20180911*
Euro OTC	Colistin sulfate	1406002-01	45	not necessary
Euro OTC	Colistin sulfate	1709009	40	20170918*
Euro OTC	Colistin sulfate	1710041	40	20171129*
Euro OTC	Colistin sulfate	1908007	40	20190905*
Euro OTC	Colistin sulfate	1905016	40	20190701*
Euro OTC	Curcuma extract ...	L1408013	60	not necessary
Euro OTC	Cyanocobalamin	1604017	10	20160509*
Euro OTC	Cyproterone acet...	1903011	20	20190404*
Euro OTC	Cyproterone acet...	1810017	20	20181108*
Euro OTC	Cyproterone acet...	1903011	40	20190404*
Euro OTC	Cyproterone acet...	1611003	30	20161201*
Euro OTC	Cysteine hydroch...	1510015	60	20151112*
Euro OTC	Cysteine hydroch...	1603018	40	20160419*
Euro OTC	Cysteine hydroch...	1304019-01	40	not necessary
Euro OTC	Cysteine hydroch...	1706023	40	20170721*
Euro OTC	Cysteine hydroch...	1507006	60	20150713*
Euro OTC	Cysteine hydroch...	1807010-01	40	not necessary
Euro OTC	Cysteine hydroch...	1901012	40	20190206*
Euro OTC	Cysteine hydroch...	1701012	40	20170207*
Euro OTC	Cystine	1811005	40	20181123*
Euro OTC	Cystine	1904022	40	20190509*
Euro OTC	D-(+)-mannose	L1503016-01	60	not necessary
Euro OTC	D-(+)-mannose	L1111005-02	40	not necessary
Euro OTC	D-(+)-mannose	L1111005-02	40	1405612
Euro OTC	Dequalinium chlo...	1405012-01	60	not necessary
Euro OTC	Dequalinium chlo...	1804009	40	20180427*
Euro OTC	Dequalinium chlo...	1512032	60	20160121*
Euro OTC	Devil's claw roo...	1210003-02	40	1707472
Euro OTC	Devil's claw roo...	1210003-02	40	not necessary
Euro OTC	Devil's claw roo...	1210003-02	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Dexamethasone	1708004	26	20170927*
Euro OTC	Dexamethasone	1806023	40	20180718*
Euro OTC	Dexamethasone	1708004	40	20170927*
Euro OTC	Dexamethasone	1303007-01	45	not necessary
Euro OTC	Dexamethasone	1806023	33	20180718*
Euro OTC	Dexamethasone	1209008-02	60	pending
Euro OTC	Dexamethasone ac...	1708011	24	20171002*
Euro OTC	Diclofenac sodium	1711026	40	20171214*
Euro OTC	Diclofenac sodium	1904008	40	20190515*
Euro OTC	Dimethyl fumarate	1601013-01	40	not necessary
Euro OTC	Dithranol	1203025-02	40	not necessary
Euro OTC	EDTA calcium dis...	1901028	40	20190218*
Euro OTC	Emulsifying ceto...	1705001	40	20170619*
Euro OTC	Emulsifying ceto...	1604013	40	20160509*
Euro OTC	Emulsifying ceto...	1701031	40	20170131*
Euro OTC	Emulsifying ceto...	1711001	40	20171121*
Euro OTC	Erythromycin	1508034	57	20150925*
Euro OTC	Erythromycin	1906032	40	20190717*
Euro OTC	Escin soluble	1612012	40	20170213*
Euro OTC	Escin soluble	1807015-01	40	not necessary
Euro OTC	Estradiol benzoate	1611004	40	20161128*
Euro OTC	Estradiol benzoate	1412033	40	20180618*
Euro OTC	Estradiol benzoate	1412033	40	20180618*
Euro OTC	Estradiol benzoate	1611004	40	20161128*
Euro OTC	Estradiol hemihy...	1808032	40	20180919*
Euro OTC	Estradiol hemihy...	1601017	45	20160122*
Euro OTC	Estradiol hemihy...	1607003	40	20160725*
Euro OTC	Estradiol hemihy...	1710031	40	20171106*
Euro OTC	Estradiol hemihy...	1906013	40	20190729*
Euro OTC	Estradiol hemihy...	1705018	30	20170630*
Euro OTC	Estradiol hemihy...	1705018	40	20170630*
Euro OTC	Estradiol hemihy...	1809022	40	20181024*
Euro OTC	Estradiol hemihy...	1808032	40	20180919*
Euro OTC	Estradiol hemihy...	1905001	40	20190521*
Euro OTC	Estradiol hemihy...	1802013	40	20180308*
Euro OTC	Estradiol hemihy...	1809022	40	20181024*
Euro OTC	Estradiol hemihy...	1906013	39	20190729*
Euro OTC	Estradiol hemihy...	1802013	40	20180308*
Euro OTC	Estradiol hemihy...	1802012	40	20180308*
Euro OTC	Estradiol hemihy...	1612008-01	30	not necessary
Euro OTC	Estradiol hemihy...	1607003	40	20160725*
Euro OTC	Estradiol hemihy...	1710031	40	20171106*
Euro OTC	Estradiol hemihy...	1905001	40	20190521*
Euro OTC	Estradiol hemihy...	1612008-01	40	not necessary
Euro OTC	Estradiol hemihy...	1802012	40	20180308*
Euro OTC	Estriol	1903035	40	20190613*
Euro OTC	Estriol	1310002-01	40	not necessary
Euro OTC	Estriol	1803011	40	20180320*
Euro OTC	Estriol	1803011	40	20180320*
Euro OTC	Estriol	1606044	40	20160803*
Euro OTC	Estriol	1610035	40	20161207*
Euro OTC	Estriol	1703044	40	20170502*
Euro OTC	Estriol	1606044	40	20160803*
Euro OTC	Estriol	1610035	40	20161207*
Euro OTC	Estriol	1710022	40	20171121*
Euro OTC	Estriol	1710030	40	20171121*
Euro OTC	Estriol	1903035	40	20190613*
Euro OTC	Estriol	1710022	40	20171121*
Euro OTC	Estriol	1710030	40	20171121*
Euro OTC	Ethyl hydrogen f...	1601011	60	20160114*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Ethyl hydrogen f...	1306014-01	60	not necessary
Euro OTC	Finasteride (for...	1809019	40	20181130*
Euro OTC	Finasteride (for...	1804010	40	20180613*
Euro OTC	Finasteride (for...	1809019	31	20181130*
Euro OTC	Finasteride (for...	1704004	40	20170426*
Euro OTC	Finasteride (for...	1804010	40	20180613*
Euro OTC	Fluocinolonaceto...	2804013E-01	46	pending
Euro OTC	Fluocinolonaceto...	2804013E-01	60	pending
Euro OTC	Folic acid hydrate	1903033	40	20190430*
Euro OTC	Folic acid hydrate	1608030	40	20160905*
Euro OTC	Frankincense ext...	1602008	60	20160215*
Euro OTC	Frankincense ext...	1308020-01	60	8958
Euro OTC	Frankincense ext...	1602008	60	20160215*
Euro OTC	Frankincense ext...	1611023-01	60	9889
Euro OTC	Fructose	1803040	40	20180410*
Euro OTC	Fructose	1710034	40	20171123*
Euro OTC	Fructose	1705002	40	20170530*
Euro OTC	Fructose	1906017	40	20190712*
Euro OTC	Fructose	1811018	40	20190205*
Euro OTC	Fructose	1607030	40	20160815*
Euro OTC	Fructose	1604024	41	20160511*
Euro OTC	Gentamicin sulfate	1604054	40	20160603*
Euro OTC	Gentamicin sulfate	1808013	40	20181008*
Euro OTC	Gentamicin sulfate	1808013	40	20181008*
Euro OTC	Gentamicin sulfate	1711025	40	20171206*
Euro OTC	Gentamicin sulfate	1506016	60	0000010321
Euro OTC	Glucose monohydr...	1908012	40	20190906*
Euro OTC	Glucose monohydr...	1210034-02	40	not necessary
Euro OTC	Glucose monohydr...	1811020	40	20181221*
Euro OTC	Glucose monohydr...	1904029	40	20190517*
Euro OTC	Glucose monohydr...	1707004	40	20170721*
Euro OTC	Glucose monohydr...	1801009	40	20180123*
Euro OTC	Glucose monohydr...	1607024	40	20160809*
Euro OTC	Glutamic acid	L1211001-01	40	not necessary
Euro OTC	Glutamic acid	1907031	40	20190815*
Euro OTC	Glutamine	1802010	40	20180215*
Euro OTC	Glutamine	1806001	40	20180614*
Euro OTC	Glutamine	1901009	40	20190215*
Euro OTC	Glycine	1906006	40	20190712*
Euro OTC	Glycine	1704015	40	20170612*
Euro OTC	Glycine	1605035	40	20160609*
Euro OTC	Griffonia seed e...	1810009	40	20181025*
Euro OTC	Guaifenesin	1701011	40	20170209*
Euro OTC	Guaifenesin	1706006	40	20170720*
Euro OTC	Guaifenesin	1901047	40	20190225*
Euro OTC	Guaifenesin	1906008	40	20190712*
Euro OTC	Guaifenesin	1804004	40	20180417*
Euro OTC	Guaifenesin	1710033	40	20171120*
Euro OTC	Guaifenesin	1809011	40	20181024*
Euro OTC	Guaifenesin	1802001	40	20180226*
Euro OTC	Heavy magnesium ...	1709011-02	40	not necessary
Euro OTC	Histamine dihydr...	1905003	40	20190521*
Euro OTC	Histamine dihydr...	1812018	40	20190129*
Euro OTC	Hyaluronic acid ...	L1302028-01	60	not necessary
Euro OTC	Hyaluronic acid ...	L1302028-01	40	1403532
Euro OTC	Hydrocortisone	1802009	40	20180302*
Euro OTC	Hydrocortisone	1611016	40	20161220*
Euro OTC	Hydrocortisone	1906033	40	20190717*
Euro OTC	Hydrocortisone	1802009	40	20180302*
Euro OTC	Hydrocortisone	1906033	40	20190717*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Hydrocortisone a...	1706021	30	20170726*
Euro OTC	Hydrocortisone a...	1805010	21	20180611*
Euro OTC	Hydrocortisone a...	1907011	40	20190813*
Euro OTC	Hydrocortisone a...	1907011	40	20190813*
Euro OTC	Hydrocortisone a...	1604052	30	20160525*
Euro OTC	Hydrocortisone a...	1805010	40	20180611*
Euro OTC	Hydrocortisone a...	1706021	40	20170726*
Euro OTC	Hydroxyethylcell...	1606030	40	20160629*
Euro OTC	Hydroxyethylcell...	1610039	40	20161128*
Euro OTC	Hydroxyethylcell...	1711012	40	20171127*
Euro OTC	Hydroxyethylcell...	1907013	40	20190725*
Euro OTC	Indometacin	1801040	40	20180215*
Euro OTC	Inositol	L1507021-01	40	not necessary
Euro OTC	Iron (II) glucon...	1805025	40	20180604*
Euro OTC	Lactose monohydr...	1802002	40	20180209*
Euro OTC	Lactose monohydr...	1703019	40	20170424*
Euro OTC	Lactose monohydr...	1812008	40	20190107*
Euro OTC	L-alanine	1710009	40	20171018*
Euro OTC	L-alanine	1601033-01	40	not necessary
Euro OTC	L-carnitine tart...	1708023	40	20170907*
Euro OTC	Levocarnitine	1303027-02	60	not necessary
Euro OTC	L-glutathione	1905011	40	20190529*
Euro OTC	L-glutathione	1711009	40	20171116*
Euro OTC	L-glutathione	1901052	40	20190226*
Euro OTC	Lidocaine	1905004	40	20190605*
Euro OTC	Lidocaine	1908005	40	20190827*
Euro OTC	Lidocaine	1801032	40	20180214*
Euro OTC	Lidocaine	1603027	40	20160413*
Euro OTC	Lidocaine	1902003	40	20190307*
Euro OTC	Lidocaine hydroc...	1606024	40	20160706*
Euro OTC	Lidocaine hydroc...	1808031	40	20180917*
Euro OTC	Lidocaine hydroc...	1612015	40	20170213*
Euro OTC	L-isoleucine	1711034	40	20171207*
Euro OTC	L-isoleucine	1708027	40	20170905*
Euro OTC	L-isoleucine	1901010	40	20190208*
Euro OTC	L-isoleucine	1805013	40	20180525*
Euro OTC	L-isoleucine	1803019	40	20180409*
Euro OTC	L-isoleucine	L1211010-02	60	not necessary
Euro OTC	L-isoleucine	L1211010-02	60	not necessary
Euro OTC	L-isoleucine	1611007	40	20161201*
Euro OTC	L-leucine	1708010	40	20170922*
Euro OTC	L-leucine	1901008	40	20190130*
Euro OTC	L-phenylalanine	1808022	40	20180830*
Euro OTC	L-phenylalanine	1907032	40	20190820*
Euro OTC	L-phenylalanine	1409013-02	40	not necessary
Euro OTC	L-proline	L1210046-02	60	not necessary
Euro OTC	L-proline	1801047	40	20180206*
Euro OTC	L-proline	L1210046-01	40	not necessary
Euro OTC	L-proline	1909002	40	20190913*
Euro OTC	L-proline	1812009	40	20190107*
Euro OTC	L-proline	1606010-01	40	1704529
Euro OTC	L-proline	L1210046-02	40	1403153
Euro OTC	L-proline	1606010-01	40	not necessary
Euro OTC	L-threonine	1711010	40	20171114*
Euro OTC	L-threonine	1609014	40	20160929*
Euro OTC	L-threonine	1703035	40	20170413*
Euro OTC	L-threonine	1606019	40	20160714*
Euro OTC	L-threonine	1811030	40	20190115*
Euro OTC	L-tryptophan	1806006	40	20180625*
Euro OTC	L-tryptophan	1906018	40	20190724*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	L-tryptophan	1812002	40	20190116*
Euro OTC	Lysine hydrochlo...	1603011	40	20160408*
Euro OTC	Lysine hydrochlo...	1902008	40	20190308*
Euro OTC	Lysine hydrochlo...	1907038	40	20190805*
Euro OTC	Macrogols 1500	1705036	40	20170627*
Euro OTC	Macrogols 1500	1711032	40	20171215*
Euro OTC	Macrogols 1500	1806015	40	20180705*
Euro OTC	Macrogols 1500	1705004	40	20170621*
Euro OTC	Macrogols 4000	1701007	40	20170213*
Euro OTC	Macrogols 4000	1803029	40	20180403*
Euro OTC	Macrogols 4000	1804001	40	20180412*
Euro OTC	Macrogols 4000	1707026	40	20170817*
Euro OTC	Macrogols 4000	1609018	40	20160929*
Euro OTC	Macrogols 4000	1810003	40	20181026*
Euro OTC	Macrogols 4000	1907002	40	20190718*
Euro OTC	Macrogols 4000	1705019	40	20170630*
Euro OTC	Macrogols 4000	1902009	40	20190313*
Euro OTC	Magnesium carbon...	1711011-01	40	not necessary
Euro OTC	Magnesium chlori...	1802023-01	40	not necessary
Euro OTC	Magnesium citrate	1701018	40	20170124*
Euro OTC	Magnesium citrate	1607006	40	20160725*
Euro OTC	Magnesium citrate	1907040	40	20190802*
Euro OTC	Magnesium citrate	1810026	40	20181114*
Euro OTC	Magnesium hydrog...	1806011	40	20180619*
Euro OTC	Magnesium hydrog...	1707010	40	20170822*
Euro OTC	Magnesium hydrog...	1707010	40	20170822*
Euro OTC	Magnesium hydrog...	1907020	40	20190805*
Euro OTC	Magnesium hydrog...	1701004	40	20170124*
Euro OTC	Magnesium hydrog...	1705010	40	20170613*
Euro OTC	Magnesium hydrog...	1904011	40	20190425*
Euro OTC	Magnesium hydrog...	1802026	40	20180305*
Euro OTC	Magnesium hydrog...	1811038	40	20190107*
Euro OTC	Magnesium hydrog...	1604002	40	20160419*
Euro OTC	Magnesium hydrog...	1608004	40	20160823*
Euro OTC	Magnesium hydrog...	1705009	40	20170613*
Euro OTC	Magnesium hydrog...	1710007	40	20171018*
Euro OTC	Magnesium orotat...	1802018	40	20180221*
Euro OTC	Magnesium orotat...	1906014	40	20190628*
Euro OTC	Magnesium orotat...	1608017	40	20160913*
Euro OTC	Magnesium orotat...	1808029	40	20180829*
Euro OTC	Magnesium orotat...	1712001	40	20171219*
Euro OTC	Magnesium orotat...	1909007	40	20191007*
Euro OTC	Magnesium orotat...	1902007	40	20190328*
Euro OTC	Magnesium stearate	1901033	40	20190320*
Euro OTC	Magnesium stearate	1608022	40	20160928*
Euro OTC	Magnesium trisil...	1608024-01	40	not necessary
Euro OTC	Manganese orotat...	1605029-01	40	not necessary
Euro OTC	Manganese orotat...	1605029-01	40	not necessary
Euro OTC	Mannitol	1703043	40	20170405*
Euro OTC	Mannitol	1904025	40	20190613*
Euro OTC	Mannitol	1902011	40	20190314*
Euro OTC	Mannitol	1610027	40	20161102*
Euro OTC	Mannitol	1908019	40	20190905*
Euro OTC	Mannitol	1708005	40	20170830*
Euro OTC	Mannitol	1606031	40	20160630*
Euro OTC	Mannitol	1604008	40	20160425*
Euro OTC	Mannitol	1809008	40	20180919*
Euro OTC	Mannitol	1806024	40	20180710*
Euro OTC	Mannitol	1802004	40	20180208*
Euro OTC	Melatonin	1811004	39	20181130*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Melatonin	1811004	40	20181130*
Euro OTC	Melatonin	1610040	40	20161201*
Euro OTC	Melatonin	1906031	40	20190709*
Euro OTC	Melatonin	1906030	40	20190709*
Euro OTC	Melatonin	1510041	60	20151123*
Euro OTC	Menthol	1701014	40	20170207*
Euro OTC	Menthol	1705022	40	20170628*
Euro OTC	Menthol	1802003	40	20180208*
Euro OTC	Menthol	1403022-02	15	not necessary
Euro OTC	Menthol	1905013	40	20190605*
Euro OTC	Mesalazine	1611019	40	20161128*
Euro OTC	Mesalazine	1902013	40	20190326*
Euro OTC	Methionine	1710008	40	20171016*
Euro OTC	Methionine	1609024	40	20161004*
Euro OTC	Methionine	1904012	40	20190509*
Euro OTC	Methionine	1703042	40	20170403*
Euro OTC	Methionine	1909004	40	20190913*
Euro OTC	Methoxsalen	1904009	40	20190628*
Euro OTC	Methoxsalen	130604101	55	not necessary
Euro OTC	Methoxsalen	1608019	34	20160913*
Euro OTC	Methoxsalen	1802021	40	20180308*
Euro OTC	Methoxsalen	1506013	45	20150706*
Euro OTC	Methyl-4-hydroxy...	1902025	40	20190313*
Euro OTC	Methyl-4-hydroxy...	1909009	40	20190913*
Euro OTC	Methyl-4-hydroxy...	1706005	40	20170726*
Euro OTC	Methyl-4-hydroxy...	1807009	40	20180718*
Euro OTC	Methylcobalamin	1711031	40	20171220*
Euro OTC	Metronidazole	1802007	40	20180223*
Euro OTC	Metronidazole	1706022	30	20170725*
Euro OTC	Metronidazole	1612013	40	20170117*
Euro OTC	Metronidazole	1901042	40	20190222*
Euro OTC	Metronidazole	1701039	40	20170307*
Euro OTC	Metronidazole	1906012	40	20190709*
Euro OTC	Metronidazole	1609028	40	20161014*
Euro OTC	Metronidazole	1807039	40	20180824*
Euro OTC	Metronidazole	1707008	30	20170725*
Euro OTC	Miconazole nitrate	1307025-02	37	not necessary
Euro OTC	Miconazole nitrate	1604004	40	20160419*
Euro OTC	Miconazole nitrate	1904032	40	20190606*
Euro OTC	Miconazole nitrate	1811015	40	20181130*
Euro OTC	Miconazole nitrate	1607004	40	20160722*
Euro OTC	Miconazole nitrate	1803028	40	20180327*
Euro OTC	Miconazole nitrate	1710020	40	20171106*
Euro OTC	Miconazole nitrate	1607004	30	20160722*
Euro OTC	Miconazole nitrate	1611024	37	20161208*
Euro OTC	Miconazole nitrate	1611024	40	20161208*
Euro OTC	Miconazole nitrate	1705039	40	20170622*
Euro OTC	Miconazole nitrate	1604004	36	20160419*
Euro OTC	Miconazole nitrate	1509014	53	20150930*
Euro OTC	Minoxidil	1602013	60	20160225*
Euro OTC	Minoxidil	1807018-01	40	not necessary
Euro OTC	Minoxidil	1509035	60	20151016*
Euro OTC	Minoxidil	1710043	40	20171122*
Euro OTC	Minoxidil	1705008	40	20170621*
Euro OTC	Minoxidil	1907025	40	20190813*
Euro OTC	Minoxidil	1411006-01	60	not necessary
Euro OTC	Minoxidil	1803004	40	20180315*
Euro OTC	Minoxidil	1607007	40	20160803*
Euro OTC	Minoxidil	1901005	40	20190128*
Euro OTC	Minoxidil	1612011	40	20170130*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Mometasone furoa...	1606021	40	20160708*
Euro OTC	Mometasone furoa...	1606021	33	20160708*
Euro OTC	Mometasone furoa...	1906035	40	20190718*
Euro OTC	Mometasone furoa...	1810008	20	20181106*
Euro OTC	Mometasone furoa...	1810008	40	20181106*
Euro OTC	Mometasone furoa...	1610009	20	20170321*
Euro OTC	Mometasone furoa...	1711027	27	20171220*
Euro OTC	Mometasone furoa...	1711027	40	20171220*
Euro OTC	Naphazoline hydr...	1903006	40	20190416*
Euro OTC	Naphazoline hydr...	1607041	40	20160819*
Euro OTC	Naphazoline hydr...	1704005	40	20170512*
Euro OTC	Naphazoline nitr...	1512033-01	40	not necessary
Euro OTC	Neomycin sulfate	1102043-02	40	not necessary
Euro OTC	Neomycin sulfate	1802005	40	20180226*
Euro OTC	Neomycin sulfate	1903022	40	20190423*
Euro OTC	Nicotinamide	1707001	40	20170711*
Euro OTC	Nicotinamide	1812019	40	20190128*
Euro OTC	Nicotinamide	1801039	40	20180205*
Euro OTC	Nicotinic acid	1801021	40	20180207*
Euro OTC	Nystatin	1612001-03	40	not necessary
Euro OTC	Ornithine aspart...	1805026	40	20180614*
Euro OTC	Ornithine aspart...	1907037	40	20190903*
Euro OTC	Ornithine aspart...	1711015	40	20171129*
Euro OTC	Ornithine hydroc...	1906015	40	20190710*
Euro OTC	Ornithine hydroc...	1704006	40	20170425*
Euro OTC	Ornithine hydroc...	1612009-01	40	not necessary
Euro OTC	Ornithine hydroc...	1803018	40	20180322*
Euro OTC	Ornithine hydroc...	1612009	40	20170113*
Euro OTC	Oxytetracycline ...	1801033	40	20180215*
Euro OTC	Pheniramine male...	1109013-02	40	not necessary
Euro OTC	Pheniramine male...	1701021	40	20170509*
Euro OTC	Phenylbutazon	1804018	40	20180504*
Euro OTC	Phenylbutazon	1604033	40	20160510*
Euro OTC	Pilocarpine hydr...	1806012	40	20180702*
Euro OTC	Pilocarpine hydr...	1806012	40	20180702*
Euro OTC	Polymyxin B sulf...	1802006	40	20180221*
Euro OTC	Polymyxin B sulf...	1807004	40	20180725*
Euro OTC	Polymyxin B sulf...	1906002	40	20190711*
Euro OTC	Polymyxin B sulf...	1701010	40	20170216*
Euro OTC	Polymyxin B sulf...	1811022	40	20190125*
Euro OTC	Polymyxin B sulf...	1202039-02	40	pending
Euro OTC	Polymyxin B sulf...	1707013	40	20170823*
Euro OTC	Potassium citrate	1603040	40	20160413*
Euro OTC	Potassium citrate	1907016	40	20190718*
Euro OTC	Potassium citrate	1606041	40	20160714*
Euro OTC	Potassium citrate	1705024	40	20170628*
Euro OTC	Potassium citrate	1902027	40	20190313*
Euro OTC	Potassium citrate	1807034	40	20180802*
Euro OTC	Potassium citrate	1803003	40	20180306*
Euro OTC	Potassium hydroc...	1805027-01	40	not necessary
Euro OTC	Potassium orotate	1203038V-02	40	9757
Euro OTC	Potassium orotate	1203038V-02	40	not necessary
Euro OTC	Potassium orotate	1203038V-02	60	not necessary
Euro OTC	Potassium perman...	1708016-02	40	not necessary
Euro OTC	Potassium sorbate	1710015	40	20171019*
Euro OTC	Potassium sorbate	1805029	40	20180611*
Euro OTC	Potassium sorbate	1703034	40	20170418*
Euro OTC	Potassium sorbate	1901020	40	20190130*
Euro OTC	Povidone iodine	1606026	40	20160708*
Euro OTC	Povidone iodine	1110013-02	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Povidone iodine	1704007	40	20170526*
Euro OTC	Prasterone	1506011	60	20150629*
Euro OTC	Prasterone	1602018	40	20160301*
Euro OTC	Prasterone	1809009	40	20181010*
Euro OTC	Prasterone	1809010	40	20181010*
Euro OTC	Prasterone	1702025	40	20170307*
Euro OTC	Prasterone	1702025	40	20170307*
Euro OTC	Prasterone	1901030	40	20190207*
Euro OTC	Prasterone	1901031	40	20190207*
Euro OTC	Prasterone	1805018	40	20180727*
Euro OTC	Prednicarbate, m...	1607012	40	20160812*
Euro OTC	Prednicarbate, m...	1805005	40	20180611*
Euro OTC	Prednicarbate, m...	1609022	27	20161024*
Euro OTC	Prednicarbate, m...	1809015	40	20181018*
Euro OTC	Prednicarbate, m...	1707022	19	20170815*
Euro OTC	Prednicarbate, m...	1605001	25	20160523*
Euro OTC	Prednicarbate, m...	1903023	40	20190503*
Euro OTC	Prednicarbate, m...	1809015	20	20181018*
Euro OTC	Prednicarbate, m...	1609022	40	20161024*
Euro OTC	Prednicarbate, m...	1903023	32	20190503*
Euro OTC	Prednicarbate, m...	1801013	20	20180126*
Euro OTC	Prednicarbate, m...	1607012	28	20160812*
Euro OTC	Prednicarbate, m...	1507011	15	20150730*
Euro OTC	Prednicarbate, m...	1801013	40	20180126*
Euro OTC	Prednicarbate, m...	1707022	40	20170815*
Euro OTC	Prednicarbate, m...	1805005	20	20180611*
Euro OTC	Prednicarbate, m...	1605001	30	20160523*
Euro OTC	Prednisolone ace...	1903031	40	20190430*
Euro OTC	Prednisolone ace...	1805009	40	20180611*
Euro OTC	Prednisolone ace...	1607010	40	20160803*
Euro OTC	Prednisolone ace...	1701041	24	20170322*
Euro OTC	Prednisolone ace...	1607010	40	20160803*
Euro OTC	Prednisolone ace...	1805009	40	20180611*
Euro OTC	Prednisolone ace...	1701041	40	20170322*
Euro OTC	Prednisolone ace...	1903031	40	20190430*
Euro OTC	Prednisolone, mi...	1702020	21	20170604*
Euro OTC	Prednisolone, mi...	1905005	40	20190613*
Euro OTC	Prednisolone, mi...	1702020	40	20170604*
Euro OTC	Prednisolone, mi...	1905005	40	20190613*
Euro OTC	Prednisone	1702008	40	20170315*
Euro OTC	Prednisone	1702008	29	20170315*
Euro OTC	Pregnenolone	1709017	40	20171010*
Euro OTC	Pregnenolone	1606028	40	20160629*
Euro OTC	Progesterone, mi...	1904033	40	20190527*
Euro OTC	Progesterone, mi...	1609016	30	20161024*
Euro OTC	Progesterone, mi...	1609001	30	20161024*
Euro OTC	Progesterone, mi...	1802030	40	20180316*
Euro OTC	Progesterone, mi...	1905022	40	20190626*
Euro OTC	Progesterone, mi...	1806020	40	20180716*
Euro OTC	Progesterone, mi...	1605013	30	20160523*
Euro OTC	Progesterone, mi...	1603012	29	20160406*
Euro OTC	Progesterone, mi...	1702026	20	20170307*
Euro OTC	Progesterone, mi...	1806020	35	20180716*
Euro OTC	Progesterone, mi...	1801036	40	20180208*
Euro OTC	Progesterone, mi...	1905023	40	20190626*
Euro OTC	Progesterone, mi...	1606033	30	20160701*
Euro OTC	Progesterone, mi...	1506033	20	0000010482
Euro OTC	Progesterone, mi...	1802030	32	20180316*
Euro OTC	Progesterone, mi...	1605032	38	20160615*
Euro OTC	Progesterone, mi...	1701038	20	20170307*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Progesterone, mi...	1603004	31	20160406*
Euro OTC	Progesterone, mi...	1905023	40	20190626*
Euro OTC	Progesterone, mi...	1905022	40	20190626*
Euro OTC	Progesterone, mi...	1511008	27	20151119*
Euro OTC	Progesterone, mi...	1606033	40	20160701*
Euro OTC	Progesterone, mi...	1705028	25	20170725*
Euro OTC	Progesterone, mi...	1603012-01	40	20181019*
Euro OTC	Progesterone, mi...	1707012	28	20170725*
Euro OTC	Progesterone, mi...	1811033	40	20181207*
Euro OTC	Progesterone, mi...	1702026	40	20170307*
Euro OTC	Progesterone, mi...	1801036	20	20180208*
Euro OTC	Progesterone, mi...	1809017	40	20181019*
Euro OTC	Progesterone, mi...	1705028	40	20170725*
Euro OTC	Progesterone, mi...	1904033	39	20190527*
Euro OTC	Progesterone, mi...	1707012	40	20170725*
Euro OTC	Progesterone, mi...	1709012	40	20171019*
Euro OTC	Progesterone, mi...	1809018	30	20181019*
Euro OTC	Progesterone, mi...	1809017	40	20181019*
Euro OTC	Progesterone, mi...	1811033	30	20181207*
Euro OTC	Progesterone, mi...	1603012-01	20	20181019*
Euro OTC	Progesterone, mi...	1609016	40	20161024*
Euro OTC	Progesterone, mi...	1603012	40	20160406*
Euro OTC	Progesterone, mi...	1605013	40	20160523*
Euro OTC	Progesterone, mi...	1601001	44	20160128*
Euro OTC	Progesterone, mi...	1506034	59	0000010483
Euro OTC	Progesterone, mi...	1509008	53	20150914*
Euro OTC	Progesterone, mi...	1709012	40	20171019*
Euro OTC	Progesterone, mi...	1901022	40	20190311*
Euro OTC	Progesterone, mi...	1901022	40	20190311*
Euro OTC	Progesterone, mi...	1901021	40	20190311*
Euro OTC	Progesterone, mi...	1701038	40	20170307*
Euro OTC	Progesterone, mi...	1605032	40	20160615*
Euro OTC	Progesterone, mi...	1901021	40	20190311*
Euro OTC	Progesterone, mi...	1609001	40	20161024*
Euro OTC	Progesterone, mi...	1809018	40	20181019*
Euro OTC	Pyridoxal-5-phos...	1808030	40	20181019*
Euro OTC	Pyridoxal-5-phos...	1711030	40	20171215*
Euro OTC	Pyridoxine hydro...	1608035	40	20160922*
Euro OTC	Pyridoxine hydro...	1806022-01	40	not necessary
Euro OTC	Pyridoxine hydro...	1604063	40	20160523*
Euro OTC	Pyridoxine hydro...	1710014	40	20171020*
Euro OTC	Pyridoxine hydro...	1908009	40	20190830*
Euro OTC	Ribavirin	1304027-01	45	AR-16-FG-007485-01
Euro OTC	Ribavirin	1304027-02	40	not necessary
Euro OTC	Riboflavin	1810005	40	20181023*
Euro OTC	Riboflavin	1809004	40	20181010*
Euro OTC	Riboflavin	1707017	40	20171012*
Euro OTC	Riboflavin	1801018	40	20180125*
Euro OTC	Ropivacaine hydr...	1707006	40	20171013*
Euro OTC	Ropivacaine hydr...	1903027	40	20190416*
Euro OTC	Salicylic acid	1901015	40	20190131*
Euro OTC	Salicylic acid	1707009	40	20170803*
Euro OTC	Salicylic acid	1801042	40	20180213*
Euro OTC	Salicylic acid	1605023	40	20160614*
Euro OTC	Selenium disulfide	1411011-02	40	not necessary
Euro OTC	Selenium yeast	L1401013-01	40	not necessary
Euro OTC	Serine	1809016	40	20181018*
Euro OTC	Serine	1611006	40	20161201*
Euro OTC	Sildenafil citrate	1902031	40	20190329*
Euro OTC	Sildenafil citrate	1807037	40	20180815*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Sodium acetate t...	1906023	40	20190628*
Euro OTC	Sodium acetate t...	1608032	40	20160914*
Euro OTC	Sodium acetate t...	1901026	40	20190204*
Euro OTC	Sodium ascorbate	1611015	40	20161202*
Euro OTC	Sodium ascorbate	1811002	40	20181207*
Euro OTC	Sodium ascorbate	1803005	40	20180314*
Euro OTC	Sodium ascorbate	1703012	40	20170313*
Euro OTC	Sodium ascorbate	1705014	40	20170630*
Euro OTC	Sodium ascorbate	1607044	40	20160823*
Euro OTC	Sodium ascorbate	1604001	40	20160413*
Euro OTC	Sodium ascorbate	1902021	40	20190313*
Euro OTC	Sodium ascorbate	1701054	40	20170307*
Euro OTC	Sodium ascorbate	1906004	40	20190709*
Euro OTC	Sodium ascorbate	1906005	40	20190709*
Euro OTC	Sodium ascorbate	1705013	40	20170630*
Euro OTC	Sodium ascorbate	1807043	40	20180802*
Euro OTC	Sodium benzoate	1803034	40	20180403*
Euro OTC	Sodium benzoate	1607005	40	20160729*
Euro OTC	Sodium benzoate	1707024	40	20170907*
Euro OTC	Sodium benzoate	1901043	40	20190206*
Euro OTC	Sodium benzoate	1809005	40	20180918*
Euro OTC	Sodium citrate	1205026-01	40	not necessary
Euro OTC	Sodium citrate	1604030	40	20160427*
Euro OTC	Sodium citrate	1909020	40	20191002*
Euro OTC	Sodium citrate	1710006	40	20171011*
Euro OTC	Sodium citrate	1607040	40	20160815*
Euro OTC	Sodium citrate	1901035	40	20190218*
Euro OTC	Sodium citrate	1807035	40	20180802*
Euro OTC	Sodium cromoglic...	1711024	40	20171219*
Euro OTC	Sodium cromoglic...	1112003-01	40	not necessary
Euro OTC	Sodium cromoglic...	1510018	60	20151109*
Euro OTC	Sodium cromoglic...	1403017-01	60	not necessary
Euro OTC	Sodium cromoglic...	1811016	40	20181218*
Euro OTC	Sodium dihydroge...	1202015-01	40	not necessary
Euro OTC	Sodium edetate	1901032	40	20190218*
Euro OTC	Sodium edetate	1704011	40	20170512*
Euro OTC	Sodium edetate	1806007	40	20180619*
Euro OTC	Sodium edetate	1609023-01	40	not necessary
Euro OTC	Sodium edetate	1807025	40	20180725*
Euro OTC	Sodium hydrogen ...	1807021-01	40	not necessary
Euro OTC	Sodium molybdate...	L1305049-01	60	not necessary
Euro OTC	Sodium molybdate...	L1305049-01	45	not necessary
Euro OTC	Sodium sulfate d...	1708012-02	40	not necessary
Euro OTC	Sodium sulfate, ...	1801046-01	40	not necessary
Euro OTC	Sodium tetraborate	1607036	40	20160816*
Euro OTC	Sodium tetraborate	1603043	40	20160419*
Euro OTC	Sodium tetraborate	1804014	40	20180426*
Euro OTC	Sodium tetraborate	1701046	40	20170131*
Euro OTC	Sorbitol	1810010	40	20181126*
Euro OTC	Spironolactone	1907001	40	20190717*
Euro OTC	Spironolactone	1711013	40	20171129*
Euro OTC	Spironolactone	1809020	40	20181015*
Euro OTC	Spironolactone	1612002	40	20170109*
Euro OTC	Sulfacetamide	1907026	40	20190814*
Euro OTC	Sulfacetamide	1702006	40	20170307*
Euro OTC	Sulfacetamide	1808035	40	20180830*
Euro OTC	Tacrolimus monoh...	1904015	40	20190617*
Euro OTC	Tacrolimus monoh...	1904015	40	20190617*
Euro OTC	Taurine	1711006	40	20171120*
Euro OTC	Taurine	1810006	40	20181116*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Testosterone	1805006	40	20180515*
Euro OTC	Testosterone	1811027	40	20190122*
Euro OTC	Testosterone	1808033	40	20180829*
Euro OTC	Testosterone	1808033	40	20180829*
Euro OTC	Testosterone	1607013	40	20160804*
Euro OTC	Testosterone	1607013	40	20160804*
Euro OTC	Testosterone	1702017	21	20170301*
Euro OTC	Testosterone	1805021	40	20180614*
Euro OTC	Testosterone	1805021	40	20180614*
Euro OTC	Testosterone	1805006	40	20180515*
Euro OTC	Testosterone	1802015	40	20180302*
Euro OTC	Testosterone	1511005-01	32	not necessary
Euro OTC	Testosterone	1511005-01	60	not necessary
Euro OTC	Testosterone	1811027	26	20190122*
Euro OTC	Testosterone	1802015	40	20180302*
Euro OTC	Testosterone	1702017	40	20170301*
Euro OTC	Testosterone pro...	1710027	40	20171114*
Euro OTC	Testosterone pro...	1811040	40	20190117*
Euro OTC	Testosterone pro...	1710027	40	20171114*
Euro OTC	Testosterone pro...	1811040	40	20190117*
Euro OTC	Tetracaine base	1903004	40	20190326*
Euro OTC	Tetracaine base	1609010	40	20161107*
Euro OTC	Tetracaine base	1805019	40	20180525*
Euro OTC	Tetracaine base	1811021	40	20181217*
Euro OTC	Tetracaine base	1708024	40	20171017*
Euro OTC	Tetracaine base	1701035	40	20170330*
Euro OTC	Tetracaine base	1606009	40	20160627*
Euro OTC	Tetracaine base	1803002	40	20180314*
Euro OTC	Tetracaine hydro...	1907007	40	20190725*
Euro OTC	Tetracaine hydro...	1607031	40	20160810*
Euro OTC	Tetracaine hydro...	1310006	40	20180116*
Euro OTC	Tetracaine hydro...	1812006	40	20190205*
Euro OTC	Tetracaine hydro...	1706018	40	20170815*
Euro OTC	Tetracaine hydro...	1712012	40	20180112*
Euro OTC	Tetracaine hydro...	1902015	40	20190326*
Euro OTC	Tetracaine hydro...	1803001	40	20180309*
Euro OTC	Tetracycline hyd...	1811031	40	20190204*
Euro OTC	Tetracycline hyd...	1706020	30	20170904*
Euro OTC	Tetracycline hyd...	1603026	40	20160607*
Euro OTC	Thiamine chlorid...	1712013	40	20180115*
Euro OTC	Thiamine chlorid...	1807001	40	20180716*
Euro OTC	Thiamine chlorid...	1705034	40	20170627*
Euro OTC	Thymol	1710037	40	20171108*
Euro OTC	Thymol	1808014	40	20180827*
Euro OTC	Titanium dioxide	1708015-01	40	not necessary
Euro OTC	Tobramycin sulfate	1710029	40	20171205*
Euro OTC	Tobramycin sulfate	1510035	60	20151208*
Euro OTC	Tobramycin sulfate	1505015-01	60	8963
Euro OTC	Tobramycin sulfate	1510035-01	60	9202
Euro OTC	Tretinoin	1612010	40	20170104*
Euro OTC	Tretinoin	1804005	40	20180423*
Euro OTC	Tretinoin	1901001	40	20190206*
Euro OTC	Tretinoin	1604059	40	20160603*
Euro OTC	Tretinoin	1612010	40	20170104*
Euro OTC	Tretinoin	1804005	40	20180423*
Euro OTC	Tretinoin	1901001	40	20190206*
Euro OTC	Triamcinolone ac...	1801020	20	20180130*
Euro OTC	Triamcinolone ac...	1810015	20	20181130*
Euro OTC	Triamcinolone ac...	1706013	20	20170720*
Euro OTC	Triamcinolone ac...	1710002	20	20171012*

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Euro OTC	Triamcinolone ac...	1508008	30	20150902*
Euro OTC	Triamcinolone ac...	1706013	40	20170720*
Euro OTC	Triamcinolone ac...	1805001-01	40	not necessary
Euro OTC	Triamcinolone ac...	1903015	40	20190410*
Euro OTC	Triamcinolone ac...	1612014	29	20170130*
Euro OTC	Triamcinolone ac...	1903015	20	20190410*
Euro OTC	Triamcinolone ac...	1810015	40	20181130*
Euro OTC	Triamcinolone ac...	1805001-01	20	not necessary
Euro OTC	Triamcinolone ac...	1801020	40	20180130*
Euro OTC	Triamcinolone ac...	1703037	40	20170526*
Euro OTC	Triamcinolone ac...	1710002	40	20171012*
Euro OTC	Triamcinolone ac...	1703037	20	20170526*
Euro OTC	Triamcinolone ac...	1612014	40	20170130*
Euro OTC	Tricalcium phosp...	1206002-01	40	1402298
Euro OTC	Triclosan	1801049	40	20180207*
Euro OTC	Triclosan	1901006	40	20190125*
Euro OTC	Triclosan	1607002	40	20160713*
Euro OTC	Triclosan	1705011	40	20170626*
Euro OTC	Triclosan	1708025	40	20170830*
Euro OTC	Tyrosine	L1401037-01	60	not necessary
Euro OTC	Ubidecarenone	1804020	40	20180529*
Euro OTC	Ubidecarenone	1810012	40	20181102*
Euro OTC	Ubidecarenone	1902030	40	20190402*
Euro OTC	Ubidecarenone	1802025	40	20180326*
Euro OTC	Urea	1702023	40	20170314*
Euro OTC	Urea	1908016	40	20190927*
Euro OTC	Urea	1702023	40	20170314*
Euro OTC	Urea	1604044	40	20160512*
Euro OTC	Urea	1707015	40	20170731*
Euro OTC	Urea	1804007	40	20180417*
Euro OTC	Urea	1811012	40	20181120*
Euro OTC	Urea	1804007	40	20180417*
Euro OTC	Urea	1811012	40	20181120*
Euro OTC	Valine	1805014	40	20180525*
Euro OTC	Valine	1705017	40	20170705*
Euro OTC	Valine	1710003	40	20171016*
Euro OTC	Valine	1609009	40	20160921*
Euro OTC	Valine	1803017	40	20180409*
Euro OTC	Valine	1712005	40	20171219*
Euro OTC	Valine	1809027	40	20181015*
Euro OTC	Vancomycin hydro...	1610038	40	20161208*
Euro OTC	Vancomycin hydro...	1604037	40	20160520*
Euro OTC	Vancomycin hydro...	1902014	40	20190321*
Euro OTC	Vancomycin hydro...	1604037	40	20160520*
Euro OTC	Vancomycin hydro...	1407047-01	60	AR-15-FG-008950-01
Euro OTC	Vancomycin hydro...	1506029	60	AR-15-FG-009077-01
Euro OTC	Vancomycin hydro...	1812005	40	20190118*
Euro OTC	Zinc gluconate	1701017	40	20170131*
Euro OTC	Zinc gluconate	1902017	40	20190517*
Euro OTC	Zinc gluconate	1810014	40	20181207*
Euro OTC	Zinc gluconate	1606018	40	20160711*
Euro OTC	Zinc gluconate	1711004	40	20171113*
Euro OTC	Zinc gluconate	1706009	40	20170727*
Euro OTC	Zinc orotate dihydro...	1203014-01	40	1407431
Euro OTC	Zinc orotate dihydro...	1203014-01	40	not necessary
Euro OTC	Zinc orotate dihydro...	13050005-02	40	not necessary
Fagron	(R,S)-epinephrin...	16H23-B98-325052	40	not necessary
Fagron	(R,S)-epinephrin...	12H08-N03	40	no certificate available
Fagron	(R,S)-epinephrin...	14A22-B05-291560	60	no certificate available
Fagron	5-aminolevulinic...	18C16-B05-347714	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	5-fluorouracil	13F19-N03	60	1608169
Fagron	5-fluorouracil	13F19-N03	40	1402308
Fagron	5-fluorouracil	16D01-U05/3-327376	40	not necessary
Fagron	8-quinolinol sul...	12H21-N08	40	not necessary
Fagron	Acetylcysteine	18B01-B04-357511	40	not necessary
Fagron	Aciclovir	15A15-B10-323362	40	not necessary
Fagron	Aciclovir	15A15-B10-323362	40	not necessary
Fagron	Aciclovir	18I18-B04-356027	40	not necessary
Fagron	Albendazole	18C21-B04-351586	40	not necessary
Fagron	Albendazole	14B05-B41-314321	40	not necessary
Fagron	Albendazole	16B04-B06-329519	40	not necessary
Fagron	Alfatradiol	15J27-B08-323877	40	not necessary
Fagron	Allopurinol	14E15-B05-330859	40	not necessary
Fagron	Allopurinol	14E15-B05-330859	40	not necessary
Fagron	Allopurinol	18E23-B02-357777	40	not necessary
Fagron	alpha-tocopheryl...	16H30-B05-326792	40	not necessary
Fagron	Aluminium chlori...	12I27-N01	40	not necessary
Fagron	Amifampridine	17J11-B01-340962	40	not necessary
Fagron	Amiloride hydroc...	15C25-B05-325653	40	not necessary
Fagron	Amiloride hydroc...	16K07-B02-329114	40	not necessary
Fagron	Amiloride hydroc...	13J08-B04-304547	60	not necessary
Fagron	Amiloride hydroc...	17F02-B06-335953	40	not necessary
Fagron	Amiloride hydroc...	13J08-B04-288620	60	AR-16-FG-008919-01
Fagron	Amitriptyline hy...	19C08-B02-362280	40	not necessary
Fagron	Amitriptyline hy...	17E11-B02-335098	40	not necessary
Fagron	Ascorbyl palmitate	17H04-B01-338484	40	not necessary
Fagron	Ascorbyl palmitate	17H04-B01-338484	40	not necessary
Fagron	Ascorbyl palmitate	18A10-B03-344364	40	not necessary
Fagron	Ascorbyl palmitate	19D25-B04-361254	40	not necessary
Fagron	Ascorbyl palmitate	18A10-B03-344364	40	not necessary
Fagron	Atropine sulfate	16G20-B06-328437	40	not necessary
Fagron	Bacitracin	16I12-B03-334194	40	not necessary
Fagron	Bacitracin	14E23-B03-296243	60	1503715
Fagron	Bacitracin	14A21-B09-292909	60	not necessary
Fagron	Bacitracin	15C12-B06-309185	60	not necessary
Fagron	Bacitracin	15C12-B06-309185	60	not necessary
Fagron	Bacitracin	17L11-B02-343685	40	1806068
Fagron	Bacitracin	18K14-B04-358690	40	not necessary
Fagron	Baclofen	18I28-B02-356294	40	not necessary
Fagron	Beclometasone di...	16B16-B03-342082	40	not necessary
Fagron	Beclometasone di...	16B16-B03-342082	25	not necessary
Fagron	Beclometasone di...	15F29-B09-314376	20	not necessary
Fagron	Beclometasone di...	15F29-B09-314376	40	not necessary
Fagron	Benzocaine	12E04-N06	30	not necessary
Fagron	Benzydamine hydr...	16D05-B50-332333	40	not necessary
Fagron	Benzydamine hydr...	16D05-B50-332333	40	not necessary
Fagron	Benzydamine hydr...	18H29-B02-355663	40	not necessary
Fagron	Benzydamine hydr...	17H01-B03-345050	40	not necessary
Fagron	Betaine hydrochl...	13H13-N02	40	1402684
Fagron	Betaine monohydr...	16H03-B04-330370	40	not necessary
Fagron	Betaine monohydr...	16H03-B04-330370	40	not necessary
Fagron	Betamethasone, m...	14C13-B02-296165	40	not necessary
Fagron	Betamethasone, m...	15L09-B05-324091	40	not necessary
Fagron	Betamethasone-17...	17F12-B10-34942	40	not necessary
Fagron	Bismuth gallath...	16K17-B04-335597	40	not necessary
Fagron	Bromelain	14F20-B03295501	40	pending
Fagron	Bromelain	14F20-B03295501	60	not necessary
Fagron	Budesonide, micr...	16J04-B01-326581	24	not necessary
Fagron	Budesonide, micr...	18C15-B03-355399	35	not necessary
Fagron	Bupivacaine hydr...	17E22-B08-345817	40	1905513

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Bupivacaine hydr...	17L05-B04-358735	40	not necessary
Fagron	Caffeine citrate	15A05-B03-318006	40	not necessary
Fagron	Caffeine citrate	15A05-B03-302746	60	not necessary
Fagron	Caffeine citrate	15A05-B03-318006	40	not necessary
Fagron	Caffeine citrate	12C16-N03	40	1403477
Fagron	Caffeine citrate	15A05-B03-302746	40	not necessary
Fagron	Calcium hydroxide	13E22-N07	20	not necessary
Fagron	Calcium lactate ...	14H01-B05-323830	40	not necessary
Fagron	Calcium orotate ...	16F13-B04-324144	40	not necessary
Fagron	Capsaicin	13D12-N01	40	1411392
Fagron	Captopril	17G21-B02-337989	30	not necessary
Fagron	Captopril	14B26-B01-312221	60	not necessary
Fagron	Chloramphenicol	16I12-B07-327124	40	not necessary
Fagron	Chlorhexidine di...	15J02-B08-322916	40	not necessary
Fagron	Chlorhexidine di...	16C11-B02-329881	40	not necessary
Fagron	Chlorhexidine di...	16C11-B02-322510	40	not necessary
Fagron	Chlorhexidine di...	14A09-B03-290724	40	AR-15-FG-012308-01
Fagron	Chlorhexidine di...	16C11-B02-324203	40	not necessary
Fagron	Chlorobutanol he...	14A06-B31-290310	60	1412458
Fagron	Chlortetracyclin...	16B10-B09-329778	30	not necessary
Fagron	Cholesterol	16G05-B03-331396	40	not necessary
Fagron	Ciclopirox olamine	16G26-B02-324884	40	not necessary
Fagron	Ciclopirox olamine	14E21-M04	60	pending
Fagron	Ciclopirox olamine	16G26-B02-324884	40	not necessary
Fagron	Ciclopirox olamine	14E21-M04	60	pending
Fagron	Ciclosporin	16C09-B02-325955	30	not necessary
Fagron	Ciclosporin	13D26-N07	18	1508429
Fagron	Cinchocaine hydr...	12F25-N10	40	1401148
Fagron	Cinchocaine hydr...	13G24-N05	40	not necessary
Fagron	Cinnarizine	11L22-N13	40	not necessary
Fagron	Cinnarizine	13H01-N07	40	1402306
Fagron	Cinnarizine	14D09-B01-294702	60	AR-16-FG-008930-01
Fagron	Ciprofloxacin hy...	16B17-B05-341882	40	not necessary
Fagron	Ciprofloxacin hy...	15D07-B01-312864	40	not necessary
Fagron	Ciprofloxacin hy...	15D07-B01-312864	40	not necessary
Fagron	Ciprofloxacin hy...	17K09-B02-345976	40	not necessary
Fagron	Citrulline	19E03-B01-362002	40	not necessary
Fagron	Clobetasol propi...	18I12-B02-357669	40	not necessary
Fagron	Clobetasol-17-pr...	13E16-N01	40	not necessary
Fagron	Clonidine hydroc...	15L11-B01-332270	40	not necessary
Fagron	Clonidine hydroc...	18B12-B04-353341	40	not necessary
Fagron	Clonidine hydroc...	14H11-B04-311163	60	AR-18-FG-006797-02
Fagron	Clonidine hydroc...	15L11-B01-332270	40	1709153
Fagron	Clonidine hydroc...	18B12-B04-353341	40	not necessary
Fagron	Clonidine hydroc...	15L11-B01-332270	40	not necessary
Fagron	Clonidine hydroc...	19B27-B02-361777	40	not necessary
Fagron	Clonidine hydroc...	15L11-B01-332270	40	not necessary
Fagron	Codeine monohydr...	12C05-N05	60	0814I-01280
Fagron	Codeine monohydr...	12C05-N05	40	0814I-01280
Fagron	Cyproheptadine h...	16K07-B01-350645	40	1906179
Fagron	Cyproterone acet...	13J16-N04	23	1602654
Fagron	Cyproterone acet...	18G25-B02-352728	20	not necessary
Fagron	Cyproterone acet...	18G25-B02-352728	40	not necessary
Fagron	Cyproterone acet...	15I11-B01-322117	33	not necessary
Fagron	Deoxycholic acid	14A02-B32-303429	40	not necessary
Fagron	Dequalinium chlo...	15B16-B04-332113	40	not necessary
Fagron	Dequalinium chlo...	11D29-N11	80	not necessary
Fagron	Dexamethasone	19A24-B01-359626	40	not necessary
Fagron	Dexamethasone	19A24-B01-359626	40	not necessary
Fagron	Dexamethasone	16F02-B05-324844	30	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Dexamethasone	19A24-B01-361445	40	not necessary
Fagron	Dexamethasone 1:...	18C13-B06-348239	40	not necessary
Fagron	Dextromethorphan...	16E23-B02-329706	40	not necessary
Fagron	Dextromethorphan...	18B09-B02-350261	40	not necessary
Fagron	D-galactose	15G13-B04-320810	40	not necessary
Fagron	D-galactose	11F30-N08	40	not necessary
Fagron	D-galactose	15G13-B04-320810	40	not necessary
Fagron	Diclofenac sodium	15D14-B04-307452	60	not necessary
Fagron	Diphenylcyclopro...	13I26-N12	40	1504598
Fagron	Diphenylcyclopro...	13I26-N12	35	not necessary
Fagron	Diphenylcyclopro...	16F02-B10-326265	40	not necessary
Fagron	Doxycycline hycl...	16I06-B02-330167	40	not necessary
Fagron	Doxycycline hycl...	18D18-B01-349567	40	not necessary
Fagron	Drofenine hydroc...	15G01-B04-330979	40	AR-17-FG-015434-02
Fagron	Drofenine hydroc...	15G01-B04-330979	40	not necessary
Fagron	Drofenine hydroc...	18A18-B99-345993	40	not necessary
Fagron	Drofenine hydroc...	15G01-B04-323028	40	not necessary
Fagron	Drofenine hydroc...	13I09-N05	40	not necessary
Fagron	Econazole nitrate	12H08-N02	40	1312442
Fagron	Econazole nitrate	12I06-B40-306925	40	AR-16-FG-011359-01
Fagron	Econazole nitrate	14G22-B07-327805	40	not necessary
Fagron	Econazole nitrate	12I06-B40-306925	40	not necessary
Fagron	Econazole nitrate	11F09-N16	40	not necessary
Fagron	Enalapril maleate	16K14-B02-336220	40	1801425
Fagron	Enalapril maleate	17J18-B03-347278	40	not necessary
Fagron	Eosin Y	13I11-N03	40	not necessary
Fagron	Eosin Y	13I11-N03	40	not necessary
Fagron	Epinephrine hydr...	18K30-B03-358514	40	not necessary
Fagron	Ergotamine tartr...	13I24-B10-288648	50	not necessary
Fagron	Ergotamine tartr...	13H02-N01	100	not necessary
Fagron	Ergotamine tartr...	13I24-B10-288650	60	not necessary
Fagron	Erythromycin	14B06-B35-290872	43	not necessary
Fagron	Estradiol hemihy...	13K22-N01	40	not necessary
Fagron	Estradiol hemihy...	16H04-B01-328738	30	not necessary
Fagron	Estriol	16I19-B01-331468	30	not necessary
Fagron	Ethinylestradiol	16B17-B01-329904	40	not necessary
Fagron	Ethinylestradiol	18E25-B01-355721	40	not necessary
Fagron	Ethinylestradiol	14E15-B03-310270	43	not necessary
Fagron	Ethinylestradiol	16B17-B01-329904	30	not necessary
Fagron	Ethyl hydrogen f...	14C12-B07-292760	60	1410436
Fagron	Ethyl hydrogen f...	14C12-B07-292760	60	not necessary
Fagron	Fampridine	16E25-B03-324097	40	not necessary
Fagron	Finasteride (for...	15J15-B02-320084	40	not necessary
Fagron	Finasteride (for...	15J15-B02-320084	40	not necessary
Fagron	Finasteride (for...	15J15-B02-320084	40	not necessary
Fagron	Finasteride (for...	15J15-B02-320084	40	not necessary
Fagron	Fluocinolonaceto...	17F14-B08-350105	40	not necessary
Fagron	Fluocinolonaceto...	14L19-B04-307101	20	not necessary
Fagron	Furosemide	16B15-B06-322860	40	not necessary
Fagron	Furosemide	16B15-B06-322860	40	not necessary
Fagron	Fusidic acid	15I10-B02-322112	40	not necessary
Fagron	Fusidic acid	17J06-B02-343675	40	not necessary
Fagron	Gentamicin sulfate	14E08-B03-294848	60	not necessary
Fagron	Gentamicin sulfate	13I19-B05-286565	60	not necessary
Fagron	Gentamicin sulfate	13I19-B05-286565	60	pending
Fagron	Glycerol monoste...	12B14-N10	40	pending
Fagron	Glycerol monoste...	11I15-N09	80	not necessary
Fagron	Glycopyrronium b...	17D28-B04-341606	40	not necessary
Fagron	Gramicidin	15H18-B01-316269	40	not necessary
Fagron	Gramicidin	14C07-B06-292886	52	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Gramicidin	14C07-B06-292886	40	not necessary
Fagron	Gramicidin	13J15-N02	60	AR-15-FG-000813-01
Fagron	Griseofulvin	17L05-B05-345995	40	not necessary
Fagron	Griseofulvin	18L04_B01-357660	40	not necessary
Fagron	Hard fat W35	15J08-B04-324959	40	not necessary
Fagron	Hard fat W35	17K09-B03-343779	40	not necessary
Fagron	Heparin sodium	16I26-B02-332108	40	not necessary
Fagron	Histamine dihydr...	15C25-B50-331242	40	not necessary
Fagron	Histamine dihydr...	15C25-B50-331242	40	not necessary
Fagron	Histidine	17F15-B04-343276	40	not necessary
Fagron	Homatropine hydr...	18L11-B03-356938	40	1905511
Fagron	Hydrochlorothiaz...	16F17-B07-327489	40	not necessary
Fagron	Hydrocortisone	16A19-B06-321981	30	not necessary
Fagron	Hydrocortisone	16A19-B06-321981	31	1702302
Fagron	Hydrocortisone a...	15A30-B09-307090	24	not necessary
Fagron	Hydrocortisone a...	18H07-B03-356259	33	not necessary
Fagron	Hydrocortisone b...	16J21-B01-333229	40	not necessary
Fagron	Hydrocortisone b...	13I26-N01	60	AR-15-FG-000944-01
Fagron	Hydrocortisone b...	15J12-B03-322394	40	not necessary
Fagron	Hydrocortisone b...	14D23-B01-296716	48	not necessary
Fagron	Hydrocortisone b...	14D23-B01-296716	60	not necessary
Fagron	Hydrocortisone b...	16J21-B01-333229	36	not necessary
Fagron	Hydrocortisone b...	16F20-B05-323903	20	not necessary
Fagron	Hydrocortisone b...	15J12-B03-322394	25	not necessary
Fagron	Hydrocortisone b...	13I26-N01	60	not necessary
Fagron	Hydrocortisone b...	15J12-B03-322394	22	not necessary
Fagron	Hydrocortisone b...	16J21-B01-333229	37	not necessary
Fagron	Hydrocortisone b...	16J21-B01-333229	40	not necessary
Fagron	Hydrocortisone b...	13I26-N01	45	not necessary
Fagron	Hydrocortisone b...	16F20-B05-323903	40	not necessary
Fagron	Hyoscine butylbr...	16I20-B08-326942	40	not necessary
Fagron	Hyoscine butylbr...	15I09-B04-316806	60	not necessary
Fagron	Hyoscine butylbr...	16I20-B08-326942	40	not necessary
Fagron	Hyoscine hydrobr...	15L09-B04-334706	40	not necessary
Fagron	Hyoscine hydrobr...	15L09-B04-334706	40	not necessary
Fagron	Indigo carmine	13K12-N01	40	1601099
Fagron	Indometacin	17I08-B06-339738	40	not necessary
Fagron	Inositol	18A03-B04-344527	40	not necessary
Fagron	Ipratropium brom...	16K07-B05-335330	40	not necessary
Fagron	Ipratropium brom...	18B19-B05-346997	40	not necessary
Fagron	Ipratropium brom...	15C04-B02-321674	40	not necessary
Fagron	Isosorbide dinit...	15I10-B01-315217	40	not necessary
Fagron	Ketamine hydroch...	15E07-B03-322164	40	not necessary
Fagron	Ketamine hydroch...	16E18-B03-339897	40	not necessary
Fagron	Ketoconazole	17H14-B01-339422	40	not necessary
Fagron	Levocarnitine	17K23-B01-345633	40	not necessary
Fagron	Levocarnitine	16F20-B01-326616	40	not necessary
Fagron	Levothyroxin sod...	17J10-B05-348764	30	not necessary
Fagron	Levothyroxin sod...	16G25-B99-329080	24	not necessary
Fagron	Levothyroxin sod...	18G27-B02-352032	32	not necessary
Fagron	Levothyroxin sod...	16G25-B99-329080	40	not necessary
Fagron	Levothyroxin sod...	17J10-B05-348764	40	not necessary
Fagron	Lithium orotate ...	18A19-B12-347442	40	not necessary
Fagron	Lithium orotate ...	16D15-B08-321402	40	not necessary
Fagron	L-malic acid	17E23-B03-335458	40	not necessary
Fagron	L-tryptophan	16E11-B03-322919	40	not necessary
Fagron	Lysine hydrochlo...	12B14-N08	40	pending
Fagron	Macrogols 6000 p...	14J17-B01-309247	60	not necessary
Fagron	Macrogols 6000 p...	14J17-B01-309247	60	AR-16-FG-011486-01
Fagron	Macrogols 6000 p...	14C11-B01-296014	60	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Mafenide acetate	15J20-B05-342526	40	not necessary
Fagron	Magnesium aspart...	18D23-B01-350069	40	not necessary
Fagron	Magnesium aspart...	15I24-B89-313231	60	not necessary
Fagron	Magnesium aspart...	17C21-B05-333111	40	not necessary
Fagron	Magnesium aspart...	15I24-B89-313231	60	AR-16-FG-008660-03
Fagron	Mannitol	18H16-B01	40	not necessary
Fagron	Melatonin	15G06-B03-315183	40	not necessary
Fagron	Melatonin	15G06-B03-323733	40	not necessary
Fagron	Mepivacaine hydr...	15I16-B01-322906	40	not necessary
Fagron	Mesalazine	17A17-B03-336444	40	not necessary
Fagron	Methacholine chl...	16L20-B01-329966	40	not necessary
Fagron	Methacholine chl...	17H03-B02-346564	40	not necessary
Fagron	Methacholine chl...	16L20-B01-329968	40	not necessary
Fagron	Methoxsalen	15J15-B03-321321	30	not necessary
Fagron	Methylprednisolone	18A03-B05-345742	40	not necessary
Fagron	Methylprednisolone	15L15-B01-321887	30	not necessary
Fagron	Methylprednisolone	18E16-B02-B02-358166	40	not necessary
Fagron	Methylprednisolone	12D17-N04	40	not necessary
Fagron	Methylprednisolone	15L15-B01-321887	40	not necessary
Fagron	Methylprednisolone	18A03-B05-345742	40	not necessary
Fagron	Methylrosanilini...	13H28-N03	60	AR-15-FG-004742-01
Fagron	Methylrosanilini...	13H28-N03	60	not necessary
Fagron	Methylrosanilini...	13G03-N11	40	not necessary
Fagron	Metoprolol tartr...	18K26-B01-359370	40	not necessary
Fagron	Metronidazole	17C10-B04-355868	40	not necessary
Fagron	Naphazoline nitr...	16A14-B03-327343	40	not necessary
Fagron	Naproxen	17C10-B03-336472	40	not necessary
Fagron	Neomycin sulfate	14C27-B05-302517	60	not necessary
Fagron	Neomycin sulfate	14C27-B05-302516	60	not necessary
Fagron	Nicotinic acid	18B15-B07-350570	40	not necessary
Fagron	Nitrofurantoin	12C30-N03	40	not necessary
Fagron	Nitrofurantoin	12C30-N03	40	1401461
Fagron	Nitrofurantoin	13B28-N07	40	not necessary
Fagron	Nitrofurantoin	11K-21-N23	40	not necessary
Fagron	Nitrofurantoin	15J15-B06-315225	60	AR-16-FG-003763-01
Fagron	Nitrofurantoin	15J15-B06-315225	60	not necessary
Fagron	Nitrofurantoin	11I07-N11	40	not necessary
Fagron	Norethisterone a...	14H29-B03-325225	40	not necessary
Fagron	Nystatin	16L21-B01-338864	40	not necessary
Fagron	Ofloxacin	15H19-B01-325226	30	not necessary
Fagron	Ofloxacin	15H19-B01-337265	40	not necessary
Fagron	Ofloxacin	15H19-B01-337265	40	not necessary
Fagron	Ofloxacin	15H19-B01-337265	40	not necessary
Fagron	Ofloxacin	15H19-B01-337265	40	not necessary
Fagron	Ofloxacin	15H19-B01-312926	60	not necessary
Fagron	Ofloxacin	13A17-N01	60	not necessary
Fagron	Omeprazole	16I20-B03-328940	40	not necessary
Fagron	Omeprazole	18K01-B03-358185	40	not necessary
Fagron	Oxybuprocaine hy...	18F20-B02-355633	40	not necessary
Fagron	Oxybuprocaine hy...	16L05-B06-334330	40	not necessary
Fagron	Oxybuprocaine hy...	16L05-B06-334330	40	not necessary
Fagron	Oxybuprocaine hy...	13J24-B07-295015	60	AR-15-FG-014812-02
Fagron	Oxybuprocaine hy...	15A26-B07-306787	60	not necessary
Fagron	Oxybutynin hydro...	12D24-N07	40	1402305
Fagron	Oxybutynin hydro...	18A12-B03-344347	40	not necessary
Fagron	Oxybutynin hydro...	16A18-B02-321181	40	not necessary
Fagron	Oxytetracycline ...	15L16-B02-322908	40	not necessary
Fagron	Oxytetracycline ...	15L16-B02-322908	40	not necessary
Fagron	Phenylephrine hy...	17C23-B04-333878	40	not necessary
Fagron	Phenylephrine hy...	17C23-B04-333880	40	not necessary

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Phenytoin	17L22-B04-345861	40	not necessary
Fagron	Phenytoin	16C14-B06-325086	40	not necessary
Fagron	Phenytoin	16C14-B06-325086	40	not necessary
Fagron	Phenytoin	12I07-N02	40	not necessary
Fagron	Polymyxin B sulf...	13E06-N12	40	pending
Fagron	Prasterone	18C12-B01-347633	40	not necessary
Fagron	Prasterone	15G14-B07-310827	60	not necessary
Fagron	Prednicarbate, m...	14A27-B02-291256	36	AR-16-FG-011442-01
Fagron	Prednisolon dihy...	16I14-B03-328918	40	not necessary
Fagron	Prednisolone 1:1...	16F09-B04-341023	40	not necessary
Fagron	Prednisolone ace...	15I23-B04-345785	40	not necessary
Fagron	Prednisolone hyd...	15C11-B05-344778	40	not necessary
Fagron	Prednisolone hyd...	12A20-N11	40	not necessary
Fagron	Prednisolone, mi...	16E09-B07-328220	40	not necessary
Fagron	Prednisolone, mi...	1108055-04	40	pending
Fagron	Prednisolone, mi...	16E09-B07-328220	20	not necessary
Fagron	Prednisone	16J19-B02-328642	30	not necessary
Fagron	Prilocaine hydro...	14H13-B03-297591	60	AR-15-FG-016323-01
Fagron	Quinidine sulfat...	16A06-B04-324954	40	not necessary
Fagron	Quinidine sulfat...	18B01-B08-350965	40	not necessary
Fagron	Quinine hydrochl...	15A20-B03-319027	40	not necessary
Fagron	Quinine hydrochl...	16K01-B05-340171	40	not necessary
Fagron	Quinoline yellow	12E24-N01	60	not necessary
Fagron	Resocrinol	17H03-B11-338899	40	not necessary
Fagron	Riboflavin phosp...	14J31-B07-300519	45	not necessary
Fagron	Riboflavin phosp...	17H03-B13-347175	40	not necessary
Fagron	Rutoside trihydr...	13K12-B06-295483	60	not necessary
Fagron	Salbutamol sulfate	13K12-N02	40	not necessary
Fagron	Salbutamol sulfate	13K12-N02	40	not necessary
Fagron	Salbutamol sulfate	13K12-N02	40	not necessary
Fagron	Salbutamol sulfate	13K12-N02	40	not necessary
Fagron	Salbutamol sulfate	17I21-B02	40	not necessary
Fagron	Salbutamol sulfate	13K12-N02	40	1801521
Fagron	Salbutamol sulfate	13K12-N02	60	not necessary
Fagron	Salbutamol sulfate	13K12-N02	40	not necessary
Fagron	Silver protein a...	18000970	40	not necessary
Fagron	Sodium alginate	14G22-B99-297091	60	not necessary
Fagron	Sodium benzoate	15B16-B07-308570	40	not necessary
Fagron	Sodium cromoglic...	15C23-B01-307751	40	not necessary
Fagron	Sodium cromoglic...	15C23-B01-307751	40	not necessary
Fagron	Sodium hydroxide...	12E02-N15	80	not necessary
Fagron	Sorbitan stearate	11I26-N09	40	not necessary
Fagron	Sunset yellow	13I05_N04	40	not necessary
Fagron	SyrSpend [®] SF Alka	17F14-B02-335945	20	1804106
Fagron	SyrSpend [®] SF Alka	15H11-B04-311722	40	1608351
Fagron	SyrSpend [®] SF Alka	15H11-B04-311722	40	not necessary
Fagron	SyrSpend [®] SF Alka	15H11-B04-311722	40	1701054
Fagron	SyrSpend [®] SF pH4	16J10-B03-330757	40	not necessary
Fagron	SyrSpend [®] SF pH4	15H13-B03-314611	40	not necessary
Fagron	SyrSpend [®] SF pH4	15H13-B03-314611	40	not necessary
Fagron	Tartaric acid	14G07-B02-322976	40	not necessary
Fagron	Testosterone	17J26-B01-342032	30	not necessary
Fagron	Tetracaine base	14A16-B04-292915	60	not necessary
Fagron	Thiamine nitrate	15A30-B06-309069	60	AR-16-FG-005855-01
Fagron	Thiomersal	11K01-N03	40	1401417
Fagron	Tranexamic acid	16E31-B07-331591	40	not necessary
Fagron	Tranexamic acid	16E31-B07-331591	50	not necessary
Fagron	Triamcinolone ac...	15C19-B06-308572	19	not necessary
Fagron	Triamcinolone ac...	13C13-N20	60	pending
Fagron	Triamcinolone, m...	17G03-B04-343670	40	1908338

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Fagron	Triamcinolone, m...	17G03-B04-343670	40	1807456
Fagron	Triamcinolone, m...	18E09-B04-353921	40	1907730
Fagron	Triamcinolone, m...	17G03-B04-343670	40	1807456
Fagron	Triamcinolone, m...	17G03-B04-343670	40	1908338
Fagron	Trometamol	17A16-B04-334121	40	not necessary
Fagron	Trometamol	17A16-B04-334121	40	not necessary
Fagron	Tropicamide	17D25-B07-336833	40	not necessary
Fagron	Tyrosine	17B06-B03-338139	40	not necessary
Fagron	Tyrothricin	17E09-B03-341231	40	not necessary
Fagron	Tyrothricin	14A23-B05-294146	46	AR-15-FG-015005-01
Fagron	Tyrothricin	12K14-N01	60	AR-15-FG-001172-01
Fagron	Zinc orotate dihydro...	13J08-N05	40	not necessary
Galke	Blackberry leaf,...	18060	40	not necessary
Galke	Boldo leaf, cut	19855	40	not necessary
Galke	Hawthorn leaf wi...	21126	40	not necessary
Galke	Maté leaf, cut	21684	40	from supplier
Galke	Melissa leaf, cut	22218	40	from supplier
Galke	Nettle leaf, cut	21238	40	from supplier
Galke	Peppermint leaf,...	20244	40	not necessary
Galke	Peppermint leaf,...	19537	40	not necessary
Galke	Silverweed, cut	22840	40	not necessary
Galke	Thyme, rubbed	21932	40	not necessary
Gatt-Koller	Caffeine citrate	4065/11164916	40	not necessary
Gatt-Koller	Caffeine citrate	1783/05140715	60	AR-15-FG-014994-01
Hepart	Acacia gum (H1001)	16000363/2	40	from supplier
Hepart	Acerola (H1413)	18000568/5	40	from supplier
Hepart	Acerola (H1413)	18001312	40	from supplier
Hepart	Amino V complex ...	18000195	40	from supplier
Hepart	Amino V complex ...	17001819	40	from supplier
Hepart	Amino V complex ...	18000193	40	from supplier
Hepart	Amino V complex ...	18000822	40	from supplier
Hepart	Amino V complex ...	18000192	40	from supplier
Hepart	Amino V complex ...	18000823	40	from supplier
Hepart	Amino V complex ...	18000824	40	from supplier
Hepart	Amino V complex ...	18000825	40	from supplier
Hepart	Amino V complex ...	18000826	40	from supplier
Hepart	Amino V complex ...	18001324	40	from supplier
Hepart	Amino V complex ...	18001575	60	from supplier
Hepart	Amino V complex ...	18001328	40	from supplier
Hepart	Amino V complex ...	18000803	40	from supplier
Hepart	Amino V complex ...	18000194	40	from supplier
Hepart	Amino V complex ...	18001674	60	from supplier
Hepart	Amino V complex ...	15001366(B1)	60	from supplier
Hepart	Amino V complex ...	18001675	60	from supplier
Hepart	Amino V complex ...	17000382	40	from supplier
Hepart	Amino V complex ...	15001366	40	from supplier
Hepart	Amino V complex ...	15001365	40	from supplier
Hepart	Amino V complex ...	16001229	40	from supplier
Hepart	Amino V complex ...	16001268	40	from supplier
Hepart	Amino V complex ...	18001676	60	from supplier
Hepart	Amino V complex ...	15001365(MM)	60	from supplier
Hepart	Amino V complex ...	15001366(F2)	60	from supplier
Hepart	Amino V complex ...	15001366(MM)	60	from supplier
Hepart	Amino V complex ...	15001365(B1)	60	from supplier
Hepart	Amino V complex ...	15001365(F2)	60	from supplier
Hepart	Aminomix NAC (H1...	041401	40	from supplier
Hepart	Aminomix NAC (H1...	041402	60	from supplier
Hepart	Aminomix NAC (H1...	14000258	40	from supplier
Hepart	Aminomix NAC (H1...	14000259	40	from supplier
Hepart	Aminomix NAC (H1...	15000298	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Aminomix NAC (H1...	14000258	40	from supplier
Hepart	Aminomix NAC (H1...	15000295 (F2)	60	from supplier
Hepart	Aminomix NAC (H1...	15000298 (F2)	60	from supplier
Hepart	Aminomix NAC (H1...	15000297	60	from supplier
Hepart	Aminomix NAC (H1...	15000321 (B1)	60	from supplier
Hepart	Aminomix NAC (H1...	15000319	60	from supplier
Hepart	Aminomix NAC (H1...	15000319 (F2)	60	from supplier
Hepart	Aminomix NAC (H1...	14000259	40	from supplier
Hepart	Aminomix NAC (H1...	041401	40	from supplier
Hepart	Aminomix NAC (H1...	15000297 (F2)	60	from supplier
Hepart	Aminomix NAC (H1...	15000295 (B1)	60	from supplier
Hepart	Aminomix NAC (H1...	17000770	40	from supplier
Hepart	Aminomix NAC (H1...	17000769	40	from supplier
Hepart	Aminomix NAC (H1...	15000295	60	from supplier
Hepart	Aminomix NAC (H1...	17000768	40	from supplier
Hepart	Aminomix NAC (H1...	17000767	40	from supplier
Hepart	Aminomix NAC (H1...	17000766	40	from supplier
Hepart	Aminomix NAC (H1...	17000771	40	from supplier
Hepart	Aminomix NAC (H1...	15000297 (B1)	60	from supplier
Hepart	Aminomix NAC (H1...	15000321	60	from supplier
Hepart	Aminomix NAC (H1...	15000320	60	from supplier
Hepart	Antibiosis forte...	17001137	40	from supplier
Hepart	Antibiosis forte...	17001151	40	from supplier
Hepart	Antibiosis forte...	15001169	60	from supplier
Hepart	Antibiosis forte...	19000078	40	from supplier
Hepart	Antibiosis mite ...	16001073	40	from supplier
Hepart	Antibiosis mite ...	16001073 (B1)	40	from supplier
Hepart	Antibiosis mite ...	16001074	40	from supplier
Hepart	Antibiosis mite ...	16001074 (B2)	40	from supplier
Hepart	Antibiosis mite ...	16001073 (B2)	40	from supplier
Hepart	Antibiosis mite ...	16001074 (B1)	40	from supplier
Hepart	Antioxidants NAT...	15000874 (B1)	60	from supplier
Hepart	Antioxidants NAT...	15000874	60	from supplier
Hepart	Antioxidants NAT...	17000260	40	from supplier
Hepart	Antioxidants NAT...	15000878 (F2)	60	from supplier
Hepart	Antioxidants NAT...	15000878 (B1)	60	from supplier
Hepart	Antioxidants NAT...	15000879	60	from supplier
Hepart	Antioxidants NAT...	15000879 (B1)	60	from supplier
Hepart	Antioxidants NAT...	15000879 (F2)	60	from supplier
Hepart	Antioxidants NAT...	15000876	60	from supplier
Hepart	Antioxidants NAT...	15000876 (B1)	60	from supplier
Hepart	Antioxidants NAT...	15000874 (F2)	60	from supplier
Hepart	Antioxidants NAT...	15000876 (F2)	60	from supplier
Hepart	Antioxidants NAT...	17000266	40	from supplier
Hepart	Antioxidants NAT...	17000263	40	from supplier
Hepart	Antioxidants NAT...	17000264	40	from supplier
Hepart	Antioxidants NAT...	14000042	40	from supplier
Hepart	Antioxidants NAT...	17000267	40	from supplier
Hepart	Antioxidants NAT...	17000269	40	from supplier
Hepart	Antioxidants NAT...	17000268	40	from supplier
Hepart	Antioxidants NAT...	19000603	40	from supplier
Hepart	Antioxidants NAT...	14000965	40	from supplier
Hepart	Antioxidants NAT...	14000961	40	from supplier
Hepart	Antioxidants NAT...	14000962	39	from supplier
Hepart	Antioxidants NAT...	14000963	40	from supplier
Hepart	Antioxidants NAT...	14000964	40	from supplier
Hepart	Antioxidants NAT...	17000262	40	from supplier
Hepart	Antioxidants NAT...	051902	60	from supplier
Hepart	Antioxidants NAT...	15000878	60	from supplier
Hepart	Antioxidants NAT...	14000966	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Antioxidants NAT...	18001818	40	from supplier
Hepart	Antioxidants NAT...	052101	50	from supplier
Hepart	Antioxidants NAT...	18001817	40	from supplier
Hepart	Antioxidants NAT...	18001816	40	from supplier
Hepart	Antioxidants NAT...	18001815	40	from supplier
Hepart	Antioxidants NAT...	18001814	40	from supplier
Hepart	Antioxidants NAT...	052201	50	from supplier
Hepart	Antioxidants NAT...	18001812	40	from supplier
Hepart	Antioxidants NAT...	052001	60	from supplier
Hepart	Antioxidants NAT...	051901	60	from supplier
Hepart	Antioxidants NAT...	18001182	40	from supplier
Hepart	Antioxidants NAT...	18001181	40	from supplier
Hepart	Antioxidants NAT...	17001755	40	from supplier
Hepart	Antioxidants NAT...	19000449	40	from supplier
Hepart	Antioxidants NAT...	14000043	40	from supplier
Hepart	Antioxidants NAT...	17001758	40	from supplier
Hepart	Antioxidants NAT...	15000875 (B1)	60	from supplier
Hepart	Antioxidants NAT...	15000875 (F2)	60	from supplier
Hepart	Antioxidants NAT...	15000880	60	from supplier
Hepart	Antioxidants NAT...	17001757	40	from supplier
Hepart	Antioxidants NAT...	17001754/0	40	from supplier
Hepart	Antioxidants NAT...	17001756	40	from supplier
Hepart	Antioxidants NAT...	067301	40	from supplier
Hepart	Antioxidants NAT...	17000261	40	from supplier
Hepart	Antioxidants NAT...	067201	40	from supplier
Hepart	Antioxidants NAT...	15000880 (B1)	60	from supplier
Hepart	Antioxidants NAT...	15000880 (F2)	60	from supplier
Hepart	Antioxidants NAT...	15000877	60	from supplier
Hepart	Antioxidants NAT...	15000877 (B1)	60	from supplier
Hepart	Antioxidants NAT...	14000044	40	from supplier
Hepart	Antioxidants NAT...	15000875	60	from supplier
Hepart	Antioxidants NAT...	15000877 (F2)	60	from supplier
Hepart	Antioxidants NAT...	17000261	40	from supplier
Hepart	Antioxidants NAT...	14001846	60	from supplier
Hepart	Antioxidants NAT...	14001851	60	from supplier
Hepart	Antioxidants NAT...	16000378	40	from supplier
Hepart	Antioxidants NAT...	16000380	40	from supplier
Hepart	Antioxidants NAT...	16000376	40	from supplier
Hepart	Antioxidants NAT...	16000377	40	from supplier
Hepart	Antioxidants NAT...	14001848	60	from supplier
Hepart	Antioxidants NAT...	14001847	60	from supplier
Hepart	Antioxidants NAT...	14001850	60	from supplier
Hepart	Antioxidants NAT...	16000974	40	from supplier
Hepart	Antioxidants NAT...	16000965	40	from supplier
Hepart	Antioxidants NAT...	16000966	40	from supplier
Hepart	Antioxidants NAT...	16000967	40	from supplier
Hepart	Antioxidants NAT...	14001849	60	from supplier
Hepart	Antioxidants NAT...	16000973	40	from supplier
Hepart	Antioxidants NAT...	16000968	40	from supplier
Hepart	Antioxidants NAT...	16000969	40	from supplier
Hepart	Antioxidants NAT...	16000972	40	from supplier
Hepart	Antioxidants NAT...	16000970	40	from supplier
Hepart	Antioxidants NAT...	16000971	40	from supplier
Hepart	Argentum (H1116)	15001151	60	from supplier
Hepart	Argentum (H1116)	16001077	40	from supplier
Hepart	Arginine (L-) (H...	066501	40	from supplier
Hepart	Arginine (L-) (H...	052401	60	from supplier
Hepart	Arginine (L-) (H...	16000949	40	from supplier
Hepart	Arginine (L-) (H...	19000626	40	from supplier
Hepart	Arginine (L-) (H...	14000563	20	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Arginine (L-) (H...	14000564	40	from supplier
Hepart	Arginine (L-) (H...	19000689	40	from supplier
Hepart	Arginine (L-) (H...	16000949	40	from supplier
Hepart	Arginine (L-) (H...	17000877	40	from supplier
Hepart	Arginine (L-) (H...	17000878	40	from supplier
Hepart	Arginine (L-) (H...	18001443	40	from supplier
Hepart	Arginine (L-) (H...	18001388	40	from supplier
Hepart	Arginine (L-) (H...	16000950	40	from supplier
Hepart	Arginine (L-) (H...	16000950	40	from supplier
Hepart	Artichoke extrac...	18001094	40	from supplier
Hepart	Baobab (H1410)	18000501/0	40	from supplier
Hepart	Baobab (H1410)	16001158	40	from supplier
Hepart	Baobab (H1410)	17001785/0	40	from supplier
Hepart	Beetroot (H1414)	18001315	40	from supplier
Hepart	Beetroot (H1414)	18000857	40	from supplier
Hepart	Betaine hydrochl...	17000829	40	from supplier
Hepart	Betaine hydrochl...	19000533	40	from supplier
Hepart	Bioflavonoid com...	068101	40	from supplier
Hepart	Bioflavonoid com...	17000257	40	from supplier
Hepart	Bioflavonoid com...	14001675	60	from supplier
Hepart	Bioflavonoid com...	17000256	40	from supplier
Hepart	Biotin (H1007)	029601	40	from supplier
Hepart	Biotin (H1007)	16001396	40	from supplier
Hepart	Biotin (H1007)	15000735	60	from supplier
Hepart	Biotin (H1007)	029601	40	from supplier
Hepart	Boron (H1005)	17000827/1	40	from supplier
Hepart	Broccoli (H1415)	18001334	40	from supplier
Hepart	Broccoli (H1415)	19000456/0	40	from supplier
Hepart	Cabbage extract ...	16000657	40	from supplier
Hepart	Cabbage extract ...	18001647	40	from supplier
Hepart	Cabbage extract ...	15001355	60	from supplier
Hepart	Cabbage extract ...	15001474	60	from supplier
Hepart	Calcium (H1008)	18001415/0	40	from supplier
Hepart	Calcium (H1008)	18001417/0	40	from supplier
Hepart	Calcium (H1008)	18001416/0	40	from supplier
Hepart	Calcium (H1008)	18001414/0	40	from supplier
Hepart	Calcium (H1008)	18001413	40	from supplier
Hepart	Calcium (H1008)	15001677(B1)	60	from supplier
Hepart	Calcium (H1008)	18001418	40	from supplier
Hepart	Calcium (H1008)	15001677/4	60	from supplier
Hepart	Calcium (H1008)	15001676(B1)	60	from supplier
Hepart	Calcium (H1008)	055101	60	from supplier
Hepart	Calcium (H1008)	15001670(B1)	60	from supplier
Hepart	Calcium (H1008)	15001665(MM)	60	from supplier
Hepart	Calcium (H1008)	15001675(MM)	60	from supplier
Hepart	Calcium (H1008)	055501	60	from supplier
Hepart	Calcium (H1008)	055401	60	from supplier
Hepart	Calcium (H1008)	055301	60	from supplier
Hepart	Calcium (H1008)	055001	60	from supplier
Hepart	Calcium (H1008)	055201	60	from supplier
Hepart	Calcium (H1008)	15001665(F2)	60	from supplier
Hepart	Calcium (H1008)	15001675(B1)	60	from supplier
Hepart	Calcium (H1008)	15001675(B2)	60	from supplier
Hepart	Calcium (H1008)	15001670(MM)	60	from supplier
Hepart	Calcium (H1008)	15001665(B1)	60	from supplier
Hepart	Calcium (H1008)	15001675/4	60	from supplier
Hepart	Calcium (H1008)	15001676(MM)	60	from supplier
Hepart	Calcium (H1008)	15001676(B2)	60	from supplier
Hepart	Calcium (H1008)	15001676/4	60	from supplier
Hepart	Calcium (H1008)	15001677(MM)	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Calcium (H1008)	055601	50	from supplier
Hepart	Calcium (H1008)	15001666 (F2)	60	from supplier
Hepart	Calcium (H1008)	15001670 (B2)	60	from supplier
Hepart	Calcium (H1008)	14000111	40	from supplier
Hepart	Calcium (H1008)	14000108	40	from supplier
Hepart	Calcium (H1008)	14000109	30	from supplier
Hepart	Calcium (H1008)	17000205	40	from supplier
Hepart	Calcium (H1008)	15001668 (B1)	60	from supplier
Hepart	Calcium (H1008)	17000204MM	40	from supplier
Hepart	Calcium (H1008)	14000106	40	from supplier
Hepart	Calcium (H1008)	15001673 (B1)	60	from supplier
Hepart	Calcium (H1008)	15001673 (B2)	60	from supplier
Hepart	Calcium (H1008)	14000107	40	from supplier
Hepart	Calcium (H1008)	15001673 (MM)	60	from supplier
Hepart	Calcium (H1008)	17000202MM	40	from supplier
Hepart	Calcium (H1008)	17000201MM	40	from supplier
Hepart	Calcium (H1008)	15001671	60	from supplier
Hepart	Calcium (H1008)	17000250	40	from supplier
Hepart	Calcium (H1008)	17000249MM	40	from supplier
Hepart	Calcium (H1008)	15001672 (MM)	60	from supplier
Hepart	Calcium (H1008)	15001672 (B1)	60	from supplier
Hepart	Calcium (H1008)	15001672 (B2)	60	from supplier
Hepart	Calcium (H1008)	15001669	60	from supplier
Hepart	Calcium (H1008)	17000203	40	from supplier
Hepart	Calcium (H1008)	15001668 (MM)	60	from supplier
Hepart	Calcium (H1008)	17000206	40	from supplier
Hepart	Calcium (H1008)	15001667 (B1)	60	from supplier
Hepart	Calcium (H1008)	15001666 (MM)	60	from supplier
Hepart	Calcium (H1008)	15001674 (B2)	60	from supplier
Hepart	Calcium (H1008)	15001667 (F2)	60	from supplier
Hepart	Calcium (H1008)	14001938	60	from supplier
Hepart	Calcium (H1008)	14001937	60	from supplier
Hepart	Calcium (H1008)	14001936	60	from supplier
Hepart	Calcium (H1008)	14001935	60	from supplier
Hepart	Calcium (H1008)	14001939	60	from supplier
Hepart	Calcium (H1008)	14001933	60	from supplier
Hepart	Calcium (H1008)	15001667 (MM)	60	from supplier
Hepart	Calcium (H1008)	14001934	60	from supplier
Hepart	Calcium (H1008)	14000110	40	from supplier
Hepart	Calcium (H1008)	15001668 (F2)	60	from supplier
Hepart	Calcium (H1008)	15001666 (B1)	60	from supplier
Hepart	Calcium (H1008)	15001674 (MM)	60	from supplier
Hepart	Calcium (H1008)	15001674 (B1)	60	from supplier
Hepart	Camu camu (H1411)	18000789/0	40	from supplier
Hepart	Camu camu (H1411)	18001526/0	40	from supplier
Hepart	Carnitine (L-) (...)	15000052	60	from supplier
Hepart	Carnitine (L-) (...)	15000053	60	from supplier
Hepart	Carnitine (L-) (...)	15000054	60	from supplier
Hepart	Carnitine (L-) (...)	15000055	60	from supplier
Hepart	Carnitine (L-) (...)	18000538	40	from supplier
Hepart	Carnitine (L-) (...)	026010	40	from supplier
Hepart	Carnitine (L-) (...)	16000453	40	from supplier
Hepart	Carnitine (L-) (...)	16000450	40	from supplier
Hepart	Carnitine (L-) (...)	16000451	40	from supplier
Hepart	Carnitine (L-) (...)	16000454	40	from supplier
Hepart	Carnitine (L-) (...)	18000418	40	from supplier
Hepart	Carnitine (L-) (...)	18000539	40	from supplier
Hepart	Carnitine (L-) (...)	130605300	40	from supplier
Hepart	Carnitine (L-) (...)	130605302	60	from supplier
Hepart	Carnosine (L-) (...)	18001109	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Carnosine (L-) (...)	074801	40	from supplier
Hepart	Carnosine (L-) (...)	074801	60	from supplier
Hepart	Carnosine (L-) (...)	045901	40	from supplier
Hepart	Carrot (H1416)	18001309	40	from supplier
Hepart	Carrot (H1416)	18001335/0	40	from supplier
Hepart	Catechin extract...	15000206/4	40	from supplier
Hepart	Catechin extract...	062901	40	from supplier
Hepart	Catechin extract...	15000206	60	from supplier
Hepart	Catechin extract...	15000206	40	from supplier
Hepart	Catechin extract...	17000409/0	40	from supplier
Hepart	Catechin extract...	062801	40	from supplier
Hepart	Catechin extract...	062701	40	from supplier
Hepart	Catechin extract...	18001420	40	from supplier
Hepart	Chaga (H1011)	16000235	60	from supplier
Hepart	Chaga (H1011)	17001560	40	from supplier
Hepart	Chaga (H1011)	16000976	40	from supplier
Hepart	Chaga (H1011)	14001688	60	from supplier
Hepart	Chaga (H1011)	14000951	40	from supplier
Hepart	Chia protein (H1...	19000342/0	40	from supplier
Hepart	Chia protein (H1...	17001729	40	from supplier
Hepart	Chondroitin sulf...	19000780	40	from supplier
Hepart	Chondroitin sulf...	14001985/2	60	from supplier
Hepart	Chondroitin sulf...	031801	60	from supplier
Hepart	Chondroitin sulf...	14001991	60	from supplier
Hepart	Chondroitin sulf...	14001985	60	from supplier
Hepart	Chondroitin sulf...	17000211	40	from supplier
Hepart	Chromium (H1012)	16000520/2	40	from supplier
Hepart	Chromium (H1012)	16000520/0	40	from supplier
Hepart	Chromium (H1012)	18002076	40	from supplier
Hepart	Chromium (H1012)	16000520/4	40	from supplier
Hepart	Chromium (H1012)	977201	15	from supplier
Hepart	Cocoa (H1424)	18001680	40	from supplier
Hepart	Cocoa (H1424)	18001645/0	40	from supplier
Hepart	Cocoa (H1424)	17000625	40	from supplier
Hepart	Coconut (H1418)	18001678	40	from supplier
Hepart	Coconut (H1418)	16001493	40	from supplier
Hepart	Coconut (H1418)	18000569	40	from supplier
Hepart	Coenzyme Q10 (H1...	16001598	40	from supplier
Hepart	Coenzyme Q10 (H1...	14001546	60	from supplier
Hepart	Coenzyme Q10 (H1...	16000402	40	from supplier
Hepart	Coenzyme Q10 (H1...	16001598/4	40	from supplier
Hepart	Coenzyme Q10 (H1...	16000402/4	40	from supplier
Hepart	Coenzyme Q10 (H1...	18000792/0	40	from supplier
Hepart	Coenzyme Q10 (H1...	16000402/0	40	from supplier
Hepart	Coenzyme Q10 (H1...	16000402/0	40	from supplier
Hepart	Coenzyme Q10 (H1...	18000806	40	from supplier
Hepart	Coenzyme Q10 (H1...	130521-300	40	from supplier
Hepart	Coenzyme Q10 (H1...	14001546	60	from supplier
Hepart	Copper (H1032)	130415	40	from supplier
Hepart	Copper (H1032)	15000211/5	40	from supplier
Hepart	Copper (H1032)	15000211	60	from supplier
Hepart	Cranberry (H1408)	18001311	40	from supplier
Hepart	Cranberry (H1408)	18000790/0	40	from supplier
Hepart	Creatine (H1031)	025003	40	from supplier
Hepart	Creatine (H1031)	025004	40	from supplier
Hepart	Creatine (H1031)	15000749	60	from supplier
Hepart	Creatine (H1031)	15000750	60	from supplier
Hepart	Creatine (H1031)	18001411	40	from supplier
Hepart	Creatine (H1031)	17001177	40	from supplier
Hepart	Curcumin and pip...	18001687	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Curcumin and pip...	18001688	40	from supplier
Hepart	Curcumin and pip...	16000369	40	from supplier
Hepart	Curcumin and pip...	14001017	40	from supplier
Hepart	Curcumin and pip...	15001589	60	from supplier
Hepart	Curcumin and pip...	17000258	40	from supplier
Hepart	Detox complex (H...	16001076	40	from supplier
Hepart	Detox complex (H...	16001075(MM)	40	from supplier
Hepart	Detox complex (H...	16001075(Nr. 5-7)	40	from supplier
Hepart	Detox complex (H...	16001075(Nr. 1-4)	40	from supplier
Hepart	Enzyme complex (...)	15001134	60	from supplier
Hepart	Flavour orange (...)	18001422/0	40	from supplier
Hepart	Flavour orange (...)	19000607/0	40	from supplier
Hepart	Flavour orange (...)	17000491/0	40	from supplier
Hepart	Flavour orange (...)	18001327	40	from supplier
Hepart	Flavour orange (...)	16000262	60	from supplier
Hepart	Flavour orange (...)	19000149	40	from supplier
Hepart	Flavour orange (...)	18000770	40	from supplier
Hepart	Folic acid (H1016)	16001098	40	from supplier
Hepart	Folic acid (H1016)	14000639	40	from supplier
Hepart	Folic acid (H1016)	052301	60	from supplier
Hepart	Folic acid (H1016)	15000125	60	from supplier
Hepart	Folic acid (H1016)	14000639	40	from supplier
Hepart	Folic acid (H1016)	14000638	40	from supplier
Hepart	Folic acid (H1016)	052301	40	from supplier
Hepart	Garcinia cambogi...	056601	40	from supplier
Hepart	Garcinia cambogi...	17000873	40	from supplier
Hepart	Garcinia cambogi...	18001849	40	from supplier
Hepart	Garcinia cambogi...	15001048(B1)	60	from supplier
Hepart	Garcinia cambogi...	15001048	60	from supplier
Hepart	Garcinia cambogi...	14001687	60	from supplier
Hepart	Garcinia cambogi...	056601	40	from supplier
Hepart	Garcinia cambogi...	15001048(F2)	60	from supplier
Hepart	Gentian extract ...	18001096	40	from supplier
Hepart	Ginger extract (...)	18001115	40	from supplier
Hepart	Ginkgo (H1024)	067901	80	from supplier
Hepart	Ginkgo (H1024)	15001038(F3)	60	from supplier
Hepart	Ginkgo (H1024)	15001039	60	from supplier
Hepart	Ginkgo (H1024)	067901	40	from supplier
Hepart	Ginkgo (H1024)	15001038	60	from supplier
Hepart	Ginkgo (H1024)	18001110	40	from supplier
Hepart	Ginkgo (H1024)	17001832	40	from supplier
Hepart	Ginkgo (H1024)	15001038(F1)	60	from supplier
Hepart	Ginkgo (H1024)	18001386	40	from supplier
Hepart	Ginseng (H1018)	018502	40	from supplier
Hepart	Ginseng (H1018)	14002047	60	from supplier
Hepart	Ginseng (H1018)	17000285	40	from supplier
Hepart	Ginseng (H1018)	17000285	40	from supplier
Hepart	Ginseng (H1018)	067801	40	from supplier
Hepart	Ginseng (H1018)	15001047	60	from supplier
Hepart	Ginseng (H1018)	067801	40	from supplier
Hepart	Ginseng (H1018)	15001047	60	from supplier
Hepart	Ginseng (H1018)	038802	40	from supplier
Hepart	Glucosamine sulf...	16000948	40	from supplier
Hepart	Glucosamine sulf...	19000775	40	from supplier
Hepart	Glucosamine sulf...	038001	40	from supplier
Hepart	Glucosamine sulf...	17001387	40	from supplier
Hepart	Glucosamine sulf...	034001	40	from supplier
Hepart	Glucosamine sulf...	17001245	40	from supplier
Hepart	Glucosamine sulf...	17000073	40	from supplier
Hepart	Glucosamine sulf...	075401	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Glucosamine sulf. ...	038001	60	from supplier
Hepart	Glutamine (L-) ...	16000942	40	from supplier
Hepart	Glutamine (L-) ...	17000901	40	from supplier
Hepart	Glutamine (L-) ...	19000650/0	40	from supplier
Hepart	Glutamine (L-) ...	17000902	40	from supplier
Hepart	Glutamine (L-) ...	19000531	40	from supplier
Hepart	Glutamine (L-) ...	15000018	60	from supplier
Hepart	Glutamine (L-) ...	061501	40	from supplier
Hepart	Glutamine (L-) ...	15000017	60	from supplier
Hepart	Glutamine (L-) ...	16000941	40	from supplier
Hepart	Glutamine (L-) ...	16000940	40	from supplier
Hepart	Glutathione (H10...	15000478	60	from supplier
Hepart	Glutathione (H10...	16000938	40	from supplier
Hepart	Glutathione (H10...	074901	40	from supplier
Hepart	Glycine (H1021)	16000933	40	from supplier
Hepart	Glycine (H1021)	16000934	40	from supplier
Hepart	Glycine (H1021)	17001658	40	from supplier
Hepart	Glycine (H1021)	17001660	40	from supplier
Hepart	Glycine (H1021)	16000238	60	from supplier
Hepart	Glycine (H1021)	14001730	60	from supplier
Hepart	Glycine (H1021)	075501	40	from supplier
Hepart	Glycine (H1021)	16000239	60	from supplier
Hepart	Glycine (H1021)	14001731	60	from supplier
Hepart	Glycine (H1021)	18000451/0	40	from supplier
Hepart	Glycine (H1021)	18000449	40	from supplier
Hepart	Glycine (H1021)	061801	40	from supplier
Hepart	Green coffee ext. ...	18000368	40	from supplier
Hepart	Hemp protein (H1...	18001330	40	from supplier
Hepart	Hemp protein (H1...	18001679	40	from supplier
Hepart	Histamine block ...	16001078	40	from supplier
Hepart	Histamine block ...	16001079	40	from supplier
Hepart	Inositol hexanic...	18001390/0	40	from supplier
Hepart	Inositol hexanic...	14000920	40	from supplier
Hepart	Inositol hexanic...	046001	60	from supplier
Hepart	Inositol hexanic...	16000750	40	from supplier
Hepart	Inositol hexanic...	15000013	60	from supplier
Hepart	Iron (H1015)	17000279	40	from supplier
Hepart	Iron (H1015)	15001258	60	from supplier
Hepart	Iron (H1015)	15001258(B1)	60	from supplier
Hepart	Iron (H1015)	14000285	60	from supplier
Hepart	Iron (H1015)	15001259	60	from supplier
Hepart	Iron (H1015)	14000285	40	from supplier
Hepart	Iron (H1015)	15001259(B1)	60	from supplier
Hepart	Iron (H1015)	15001259(F2)	60	from supplier
Hepart	Iron (H1015)	15001258(F2)	60	from supplier
Hepart	Iron (H1015)	14000285	60	from supplier
Hepart	Iron (H1015)	040901	60	from supplier
Hepart	Iron (H1015)	040901	60	from supplier
Hepart	Isoflavones (H10...	16001644/0	40	from supplier
Hepart	Isoflavones (H10...	15000205	60	from supplier
Hepart	Isoflavones (H10...	15000205/2	40	from supplier
Hepart	Isoflavones (H10...	056701	40	from supplier
Hepart	Isoleucine (L-) ...	14001772	60	from supplier
Hepart	Isoleucine (L-) ...	17000900	40	from supplier
Hepart	Isoleucine (L-) ...	15001641	60	from supplier
Hepart	Isoleucine (L-) ...	18001412	40	from supplier
Hepart	Isoleucine (L-) ...	15001640	60	from supplier
Hepart	Isoleucine (L-) ...	066301	40	from supplier
Hepart	Juniper (H1112)	16000683	40	from supplier
Hepart	Juniper (H1112)	15001004	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Konjac (H1426)	17000846	40	from supplier
Hepart	Konjac (H1426)	19000132	40	from supplier
Hepart	Konjac (H1426)	17000837	40	from supplier
Hepart	L-5-HTP (H1035)	16000441	40	from supplier
Hepart	L-5-HTP (H1035)	16000441/0	40	from supplier
Hepart	Leaky gut comple...	18000639	40	from supplier
Hepart	Leaky gut comple...	18000638	40	from supplier
Hepart	Leucine (L-) (H1...	18001305	40	from supplier
Hepart	Leucine (L-) (H1...	074301	40	from supplier
Hepart	Leucine (L-) (H1...	965101	60	from supplier
Hepart	Leucine (L-) (H1...	15000283	60	from supplier
Hepart	Lipoic acid (a-)...	14001726	60	from supplier
Hepart	Lipoic acid (a-)...	054301	58	from supplier
Hepart	Lipoic acid (a-)...	18000250	40	from supplier
Hepart	Lipoic acid (a-)...	14000497 (F2)	60	from supplier
Hepart	Lipoic acid (a-)...	14001420	60	from supplier
Hepart	Lipoic acid (a-)...	15001526 (B1.1)	60	from supplier
Hepart	Lipoic acid (a-)...	14000497	60	from supplier
Hepart	Lipoic acid (a-)...	066601	40	from supplier
Hepart	Lipoic acid (a-)...	15001526 (B1.2)	60	from supplier
Hepart	Lipoic acid (a-)...	14000497 (F1)	60	from supplier
Hepart	Lipoic acid (a-)...	15001526 (MM)	60	from supplier
Hepart	Liver kidney com...	18000635	40	from supplier
Hepart	Lupin protein co...	17001732	40	from supplier
Hepart	Lupin protein co...	18000948	40	from supplier
Hepart	Lupin protein co...	19000532	40	from supplier
Hepart	Lycopene (H1058)	17000619/0	40	from supplier
Hepart	Lycopene (H1058)	17000619	40	from supplier
Hepart	Lycopene (H1058)	053501	40	from supplier
Hepart	Lysine (L-) (H10...	065401	40	from supplier
Hepart	Lysine (L-) (H10...	028801	40	from supplier
Hepart	Lysine (L-) (H10...	16001503/0	40	from supplier
Hepart	Lysine (L-) (H10...	14000944	40	from supplier
Hepart	Lysine (L-) (H10...	18001389	40	from supplier
Hepart	Lysine (L-) (H10...	065401	40	from supplier
Hepart	Lysine (L-) (H10...	15000372	60	from supplier
Hepart	Lysine (L-) (H10...	16001395/0	40	from supplier
Hepart	Lysine (L-) (H10...	028801	40	from supplier
Hepart	Magnesium (H1060)	18000953	40	from supplier
Hepart	Magnesium (H1060)	19001559/0	40	from supplier
Hepart	Magnesium (H1060)	14001889	60	from supplier
Hepart	Magnesium (H1060)	16000272	60	from supplier
Hepart	Magnesium (H1060)	18000952	40	from supplier
Hepart	Magnesium (H1060)	18000951	40	from supplier
Hepart	Magnesium (H1060)	18000950	40	from supplier
Hepart	Magnesium (H1060)	18000949	40	from supplier
Hepart	Magnesium (H1060)	18000420	40	from supplier
Hepart	Magnesium (H1060)	19001560	40	from supplier
Hepart	Magnesium (H1060)	19001153	40	from supplier
Hepart	Magnesium (H1060)	16000269	60	from supplier
Hepart	Magnesium (H1060)	19001561	40	from supplier
Hepart	Magnesium (H1060)	069701	40	from supplier
Hepart	Magnesium (H1060)	14001895	60	from supplier
Hepart	Magnesium (H1060)	043601	40	from supplier
Hepart	Magnesium (H1060)	16000271	60	from supplier
Hepart	Magnesium (H1060)	14001893	60	from supplier
Hepart	Magnesium (H1060)	069601	40	from supplier
Hepart	Magnesium (H1060)	16001705	40	from supplier
Hepart	Magnesium (H1060)	14001891	60	from supplier
Hepart	Magnesium (H1060)	14001894	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Magnesium (H1060)	16000265	60	from supplier
Hepart	Magnesium (H1060)	16001704	40	from supplier
Hepart	Magnesium (H1060)	16000270	60	from supplier
Hepart	Magnesium (H1060)	16000264	60	from supplier
Hepart	Magnesium (H1060)	16001706	40	from supplier
Hepart	Magnesium (H1060)	16000268	60	from supplier
Hepart	Magnesium (H1060)	070001	40	from supplier
Hepart	Magnesium (H1060)	19001558	40	from supplier
Hepart	Magnesium (H1060)	069801	40	from supplier
Hepart	Magnesium (H1060)	16000266	60	from supplier
Hepart	Magnesium (H1060)	069901	40	from supplier
Hepart	Magnesium (H1060)	16000267	60	from supplier
Hepart	Magnesium (H1060)	16001703	40	from supplier
Hepart	Magnesium (H1060)	16001702	40	from supplier
Hepart	Marshmallow extr...	18001095	40	from supplier
Hepart	Melatonin (H1061)	069101	60	from supplier
Hepart	Melatonin (H1061)	069101	40	from supplier
Hepart	Methylsulfonylme...	18001387	40	from supplier
Hepart	Methylsulfonylme...	16000449/0	40	from supplier
Hepart	Methylsulfonylme...	056201	60	from supplier
Hepart	Methylsulfonylme...	17001247	40	from supplier
Hepart	Methylsulfonylme...	15001371	120	from supplier
Hepart	Methylsulfonylme...	14000921	40	from supplier
Hepart	OPC Grape seed (...)	17001805	40	from supplier
Hepart	OPC Grape seed (...)	16001439(B1)	40	from supplier
Hepart	OPC Grape seed (...)	16001439	40	from supplier
Hepart	OPC Grape seed (...)	16001439(F1)	40	from supplier
Hepart	OPC Grape seed (...)	068301	40	from supplier
Hepart	OPC Grape seed (...)	070401	40	from supplier
Hepart	OPC Grape seed (...)	070401	40	from supplier
Hepart	OPC Grape seed (...)	130208-77	40	from supplier
Hepart	Ornithine (L-) (...)	15001364(F1)	60	from supplier
Hepart	Ornithine (L-) (...)	15001364(F2)	60	from supplier
Hepart	Ornithine (L-) (...)	15000004	60	from supplier
Hepart	Ornithine (L-) (...)	16000935	40	from supplier
Hepart	Ornithine (L-) (...)	15001364	60	from supplier
Hepart	Orthocolon compl...	18000636	40	from supplier
Hepart	Orthocolon compl...	19000937	40	from supplier
Hepart	Orthocolon compl...	18000640	40	from supplier
Hepart	Orthocolon compl...	18000641	40	from supplier
Hepart	Orthovimin Neutr...	17001152	40	from supplier
Hepart	Orthovimin Neutr...	19001312	40	from supplier
Hepart	Orthovimin Neutr...	18001065	40	from supplier
Hepart	Orthovimin Neutr...	15001850	60	from supplier
Hepart	Orthovimin Neutr...	17001625	40	from supplier
Hepart	Orthovimin Neutr...	19000033	40	from supplier
Hepart	Pathogene comple...	18000637	40	from supplier
Hepart	Pea protein comp...	18001020	40	from supplier
Hepart	Pea protein comp...	18001019	40	from supplier
Hepart	Pea protein comp...	17001708	40	from supplier
Hepart	Pea protein comp...	17001707	40	from supplier
Hepart	PET basis (H1503)	15000402	60	from supplier
Hepart	PET basis (H1503)	15000402	60	from supplier
Hepart	PET coenzyme Q10...	15000401	60	from supplier
Hepart	PET coenzyme Q10...	15000401	60	from supplier
Hepart	PET immune syste...	15000405	60	from supplier
Hepart	PET immune syste...	15000405	60	from supplier
Hepart	PET joints (H1501)	15000404	60	from supplier
Hepart	PET joints (H1501)	15000404	60	from supplier
Hepart	PET skin and hai...	15000403	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	PET skin and hai...	15000403	60	from supplier
Hepart	PET well-aging (...)	15000406	60	from supplier
Hepart	PET well-aging (...)	15000401	60	from supplier
Hepart	Phenylalanine (L...	19000779	40	from supplier
Hepart	Phenylalanine (L...	17000876	40	from supplier
Hepart	Phenylalanine (L...	16001096	40	from supplier
Hepart	Phenylalanine (L...	074401	40	from supplier
Hepart	Pomegranate (H14...	18001310	40	from supplier
Hepart	Pomegranate (H14...	18000566/0	40	from supplier
Hepart	Potassium (H1030)	16001646	40	from supplier
Hepart	Potassium (H1030)	16001647	40	from supplier
Hepart	Potassium (H1030)	15000599	60	from supplier
Hepart	Potassium (H1030)	16001645	40	from supplier
Hepart	Potassium (H1030)	15000599(F2)	60	from supplier
Hepart	Potassium (H1030)	15000600	60	from supplier
Hepart	Potassium (H1030)	15000597	60	from supplier
Hepart	Potassium (H1030)	15000598	60	from supplier
Hepart	Potassium (H1030)	15000598(F2)	60	from supplier
Hepart	Potassium (H1030)	19000028	40	from supplier
Hepart	Potassium (H1030)	19000027	40	from supplier
Hepart	Potassium (H1030)	19000026	40	from supplier
Hepart	Potassium (H1030)	130328001	40	from supplier
Hepart	Potassium (H1030)	15000597(B1)	60	from supplier
Hepart	Potassium (H1030)	069201	30	from supplier
Hepart	Potassium (H1030)	039302	40	from supplier
Hepart	Potassium (H1030)	15000600(B1)	60	from supplier
Hepart	Potassium (H1030)	15000598(B1)	60	from supplier
Hepart	Proline (-L) (H1...	066401	40	from supplier
Hepart	Proline (-L) NF ...	16000951	40	from supplier
Hepart	Proline (-L) NF ...	17000952	40	from supplier
Hepart	Proline (-L) NF ...	16000951	40	from supplier
Hepart	Proline (-L) NF ...	18000403	40	from supplier
Hepart	Proline (-L) NF ...	16000951	40	from supplier
Hepart	Psyllium husks (...)	19000338/0	40	from supplier
Hepart	Psyllium husks (...)	17001728	40	from supplier
Hepart	Pumpkin seed pro...	17001731	40	from supplier
Hepart	Pumpkin seed pro...	19000340/0	40	from supplier
Hepart	Rhodiola (H1064)	15001045	60	from supplier
Hepart	Rhodiola (H1064)	070501	60	from supplier
Hepart	Rhodiola (H1064)	18000413	40	from supplier
Hepart	Rhodiola (H1064)	070501	40	from supplier
Hepart	Rhodiola (H1064)	17000278	40	from supplier
Hepart	Saffron flowers ...	18000806	40	from supplier
Hepart	Saffron style (H...	18001093	40	from supplier
Hepart	Seabuckthorn org...	16001709	40	from supplier
Hepart	Seabuckthorn org...	16001154	40	from supplier
Hepart	Seabuckthorn org...	18001314	40	from supplier
Hepart	Selenium (H1067)	18000422	40	from supplier
Hepart	Selenium (H1067)	065501	40	from supplier
Hepart	Selenium (H1067)	18001022	40	from supplier
Hepart	Selenium (H1067)	18001021	40	from supplier
Hepart	Selenium (H1067)	16001349	40	from supplier
Hepart	Selenium (H1067)	18000578	40	from supplier
Hepart	Selenium (H1067)	065601	40	from supplier
Hepart	Selenium (H1067)	16001350	40	from supplier
Hepart	Selenium (H1067)	16001351	40	from supplier
Hepart	Selenium (H1067)	15000157	60	from supplier
Hepart	Selenium (H1067)	15000159	60	from supplier
Hepart	Selenium (H1067)	15000158	60	from supplier
Hepart	Selenium (H1067)	065701	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Selenium (H1067)	065501	60	from supplier
Hepart	Selenium yeast (...)	16001450/1	40	from supplier
Hepart	Silicium (H1070)	18000448	40	from supplier
Hepart	Silicium (H1070)	18001681	40	from supplier
Hepart	Silicium (H1070)	17000326	40	from supplier
Hepart	Silicium (H1070)	17000817	40	from supplier
Hepart	Silicium (H1070)	19000534	40	from supplier
Hepart	Silymarin (H1071)	067701	40	from supplier
Hepart	Silymarin (H1071)	16001084	40	from supplier
Hepart	Silymarin (H1071)	16001084/0	40	from supplier
Hepart	Silymarin (H1071)	15001046	60	from supplier
Hepart	Stomach complex ...	16001139	40	from supplier
Hepart	Stomach complex ...	15001350	60	from supplier
Hepart	Stomach complex ...	16001138	40	from supplier
Hepart	Sunflower seed p...	19000339/0	40	from supplier
Hepart	Sunflower seed p...	17001730	40	from supplier
Hepart	Taurine (H1084)	18000419	40	from supplier
Hepart	Taurine (H1084)	16000946	40	from supplier
Hepart	Taurine (H1084)	16000946	40	from supplier
Hepart	Taurine (H1084)	14001767	60	from supplier
Hepart	Taurine (H1084)	19000847	40	from supplier
Hepart	Taurine (H1084)	16000947	40	from supplier
Hepart	Threonine (L-) (...)	18001023	40	from supplier
Hepart	Threonine (L-) (...)	14001718	40	from supplier
Hepart	Threonine (L-) (...)	14001718	60	from supplier
Hepart	Threonine (L-) (...)	17000875	40	from supplier
Hepart	TopFlora multi f...	15001135	60	from supplier
Hepart	Trace elements B...	16000796	40	from supplier
Hepart	Trace elements B...	16000795	40	from supplier
Hepart	Trace elements B...	16000797	40	from supplier
Hepart	Trace elements J...	15000525	60	from supplier
Hepart	Trace elements J...	18000456	40	from supplier
Hepart	Trace elements J...	18000457	40	from supplier
Hepart	Trace elements J...	072901	60	from supplier
Hepart	Trace elements J...	18000484	40	from supplier
Hepart	Trace elements J...	15000526	60	from supplier
Hepart	Trace elements J...	072701	40	from supplier
Hepart	Trace elements J...	17000870	40	from supplier
Hepart	Trace elements J...	072801	40	from supplier
Hepart	Trace elements J...	072901	40	from supplier
Hepart	Trace elements J...	17000871	40	from supplier
Hepart	Trace elements J...	17000869	40	from supplier
Hepart	Trace elements J...	18000455	40	from supplier
Hepart	Trace elements S...	042301	15	from supplier
Hepart	Trace elements S...	14000721	40	from supplier
Hepart	Trace elements S...	19001256	40	from supplier
Hepart	Trace elements S...	14000722	40	from supplier
Hepart	Trace elements S...	14000720	40	from supplier
Hepart	Trace elements S...	14000721	40	from supplier
Hepart	Trace elements S...	14000722	40	from supplier
Hepart	Trace elements S...	14000720	40	from supplier
Hepart	Trace elements S...	19001260	40	from supplier
Hepart	Trace elements S...	15000555	60	from supplier
Hepart	Trace elements S...	042301	15	from supplier
Hepart	Trace elements S...	15000558(B1)	60	from supplier
Hepart	Trace elements S...	15000557(F2)	60	from supplier
Hepart	Trace elements S...	15000558	60	from supplier
Hepart	Trace elements S...	15000557	60	from supplier
Hepart	Trace elements S...	15000556(F2)	60	from supplier
Hepart	Trace elements S...	15000558(F2)	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Trace elements S...	15000556(B1)	60	from supplier
Hepart	Trace elements S...	15000558	40	from supplier
Hepart	Trace elements S...	17000886	40	from supplier
Hepart	Trace elements S...	15000557(B1)	60	from supplier
Hepart	Trace elements S...	17000888	40	from supplier
Hepart	Trace elements S...	17000887	40	from supplier
Hepart	Trace elements S...	15000555(B1)	60	from supplier
Hepart	Trace elements S...	15000555(F2)	60	from supplier
Hepart	Trace elements S...	15000553	60	from supplier
Hepart	Trace elements S...	15000553(F2)	60	from supplier
Hepart	Trace elements S...	15000553(B1)	60	from supplier
Hepart	Trace elements S...	15000556	60	from supplier
Hepart	Tryptophan (L-) ...	074501	40	from supplier
Hepart	Tryptophan (L-) ...	16001394	40	from supplier
Hepart	Tryptophan (L-) ...	15000443	60	from supplier
Hepart	Tryptophan (L-) ...	15000441	60	from supplier
Hepart	Tryptophan (L-) ...	16001404	40	from supplier
Hepart	Tyrosine (L-) (H...	074601	100	from supplier
Hepart	Tyrosine (L-) (H...	15001380	60	from supplier
Hepart	Tyrosine (L-) (H...	18000412	40	from supplier
Hepart	Valine (L-) (H10...	15001436	60	from supplier
Hepart	Valine (L-) (H10...	14001947	60	from supplier
Hepart	Valine (L-) (H10...	16001380/0	40	from supplier
Hepart	Vitamin B1 (H1101)	075001	60	from supplier
Hepart	Vitamin B1 (H1101)	17001075	40	from supplier
Hepart	Vitamin B1 (H1101)	17000914	40	from supplier
Hepart	Vitamin B1 (H1101)	19000685	40	from supplier
Hepart	Vitamin B1 (H1101)	053301	60	from supplier
Hepart	Vitamin B1 (H1101)	075001	40	from supplier
Hepart	Vitamin B1 (H1101)	14001112	40	from supplier
Hepart	Vitamin B1 (H1101)	14001113	40	from supplier
Hepart	Vitamin B1 (H1101)	075001	40	from supplier
Hepart	Vitamin B12 (H11...	17000277	40	from supplier
Hepart	Vitamin B12 (H11...	17000276	40	from supplier
Hepart	Vitamin B12 (H11...	15000993	60	from supplier
Hepart	Vitamin B12 (H11...	025017	40	from supplier
Hepart	Vitamin B12 (H11...	15000993(F2)	60	from supplier
Hepart	Vitamin B12 (H11...	15000993(B1)	60	from supplier
Hepart	Vitamin B6 (H1103)	19000782	40	from supplier
Hepart	Vitamin B6 (H1103)	14001576	40	from supplier
Hepart	Vitamin B6 (H1103)	16001239	40	from supplier
Hepart	Vitamin B6 (H1103)	16001239	40	from supplier
Hepart	Vitamin B6 (H1103)	15000279	60	from supplier
Hepart	Vitamin B6 (H1103)	15000279	60	from supplier
Hepart	Vitamin C (H1104)	19001006	40	from supplier
Hepart	Vitamin C (H1104)	074201	40	from supplier
Hepart	Vitamin C (H1104)	19000134	40	from supplier
Hepart	Vitamin C (H1104)	16000199	60	from supplier
Hepart	Vitamin C (H1104)	16000198	60	from supplier
Hepart	Vitamin C (H1104)	19001005	40	from supplier
Hepart	Vitamin C (H1104)	16000197	60	from supplier
Hepart	Vitamin C (H1104)	074201	38	from supplier
Hepart	Vitamin C (H1104)	074001	40	from supplier
Hepart	Vitamin C (H1104)	19001003	40	from supplier
Hepart	Vitamin C (H1104)	17000920	40	from supplier
Hepart	Vitamin C (H1104)	14001794	60	from supplier
Hepart	Vitamin C (H1104)	074101	40	from supplier
Hepart	Vitamin C (H1104)	16000203	60	from supplier
Hepart	Vitamin C (H1104)	16000202	60	from supplier
Hepart	Vitamin C (H1104)	073901	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Vitamin C (H1104)	16000231	60	from supplier
Hepart	Vitamin C (H1104)	16000230	60	from supplier
Hepart	Vitamin C (H1104)	044602	40	from supplier
Hepart	Vitamin C (H1104)	16000229	60	from supplier
Hepart	Vitamin C (H1104)	14001797	60	from supplier
Hepart	Vitamin C (H1104)	16000201	60	from supplier
Hepart	Vitamin C (H1104)	14001788	60	from supplier
Hepart	Vitamin C (H1104)	14001780	60	from supplier
Hepart	Vitamin C (H1104)	17000915	40	from supplier
Hepart	Vitamin C (H1104)	17000916	40	from supplier
Hepart	Vitamin C (H1104)	17000917	40	from supplier
Hepart	Vitamin C (H1104)	044502	40	from supplier
Hepart	Vitamin C (H1104)	17000918	40	from supplier
Hepart	Vitamin C (H1104)	17000919	40	from supplier
Hepart	Vitamin C (H1104)	16000200	60	from supplier
Hepart	Vitamin C (H1104)	19001004	40	from supplier
Hepart	Vitamin complex ...	16000960	40	from supplier
Hepart	Vitamin complex ...	19001423/0	40	from supplier
Hepart	Vitamin complex ...	19001424/0	40	from supplier
Hepart	Vitamin complex ...	19000433	40	from supplier
Hepart	Vitamin complex ...	19000544	40	from supplier
Hepart	Vitamin complex ...	14001590	60	from supplier
Hepart	Vitamin complex ...	14001589	60	from supplier
Hepart	Vitamin complex ...	14001591	60	from supplier
Hepart	Vitamin complex ...	16000957	40	from supplier
Hepart	Vitamin complex ...	16000958	40	from supplier
Hepart	Vitamin complex ...	16000959	40	from supplier
Hepart	Vitamin complex ...	16000961	40	from supplier
Hepart	Vitamin complex ...	049201	40	from supplier
Hepart	Vitamin complex ...	16000962	40	from supplier
Hepart	Vitamin complex ...	17000880/0	40	from supplier
Hepart	Vitamin complex ...	15001276(B1)	60	from supplier
Hepart	Vitamin complex ...	15001276	60	from supplier
Hepart	Vitamin complex ...	15001276(F2)	60	from supplier
Hepart	Vitamin complex ...	15001275(F2)	60	from supplier
Hepart	Vitamin complex ...	15001274	60	from supplier
Hepart	Vitamin complex ...	15001277(B1)	60	from supplier
Hepart	Vitamin complex ...	15001274(B1)	60	from supplier
Hepart	Vitamin complex ...	15001275(B1)	60	from supplier
Hepart	Vitamin complex ...	15001275	60	from supplier
Hepart	Vitamin complex ...	15001277(F2)	60	from supplier
Hepart	Vitamin complex ...	077101	40	from supplier
Hepart	Vitamin complex ...	15001277	60	from supplier
Hepart	Vitamin complex ...	17001753/0	40	from supplier
Hepart	Vitamin complex ...	077001	40	from supplier
Hepart	Vitamin complex ...	15001274(F2)	60	from supplier
Hepart	Vitamin complex ...	17001752/0	40	from supplier
Hepart	Vitamin complex ...	17001751/0	40	from supplier
Hepart	Vitamin complex ...	17001750/0	40	from supplier
Hepart	Vitamin complex ...	17001654	40	from supplier
Hepart	Vitamin complex ...	077201	40	from supplier
Hepart	Vitamin D3 (H1105)	15000490	60	from supplier
Hepart	Vitamin D3 (H1105)	16001212	40	from supplier
Hepart	Vitamin D3 (H1105)	17000874	40	from supplier
Hepart	Vitamin D3 (H1105)	18001951	40	from supplier
Hepart	Vitamin D3 (H1105)	16000952	40	from supplier
Hepart	Vitamin E NAT (H...	17000890	40	from supplier
Hepart	Vitamin E NAT (H...	15001559	60	from supplier
Hepart	Vitamin E NAT (H...	15001559(F2)	60	from supplier
Hepart	Vitamin E NAT (H...	071601	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Hepart	Vitamin E NAT (H...	15001561	60	from supplier
Hepart	Vitamin E NAT (H...	15001560(B1)	60	from supplier
Hepart	Vitamin E NAT (H...	15001561(B1)	60	from supplier
Hepart	Vitamin E NAT (H...	071701	35	from supplier
Hepart	Vitamin E NAT (H...	065001	40	from supplier
Hepart	Vitamin E NAT (H...	15001561(F2)	60	from supplier
Hepart	Vitamin E NAT (H...	17000892	40	from supplier
Hepart	Vitamin E NAT (H...	15001558(F2)	60	from supplier
Hepart	Vitamin E NAT (H...	15001558(B1)	60	from supplier
Hepart	Vitamin E NAT (H...	15001557(B1)	60	from supplier
Hepart	Vitamin E NAT (H...	16001246	40	from supplier
Hepart	Vitamin E NAT (H...	15001558	60	from supplier
Hepart	Vitamin E NAT (H...	15001560(F2)	60	from supplier
Hepart	Vitamin E NAT (H...	15001559(B1)	60	from supplier
Hepart	Vitamin E NAT (H...	15001560	60	from supplier
Hepart	Vitamin E NAT (H...	15001557(F2)	60	from supplier
Hepart	Vitamin E NAT (H...	071501	40	from supplier
Hepart	Vitamin E NAT (H...	15001557	60	from supplier
Hepart	Vitamin E NAT (H...	14001821	60	from supplier
Hepart	Vitamin E NAT (H...	14001822	60	from supplier
Hepart	Vitamin E NAT (H...	071501	40	from supplier
Hepart	Vitamin E NAT (H...	071501	40	from supplier
Hepart	Vitamin E NAT (H...	071501	40	from supplier
Hepart	Vitamin E NAT (H...	14001815	60	from supplier
Hepart	Vitamin E NAT (H...	18001026	40	from supplier
Hepart	Vitamin E NAT (H...	18001025	40	from supplier
Hepart	Vitamin E NAT (H...	18001024	40	from supplier
Hepart	Vitamin E NAT (H...	17000891	40	from supplier
Hepart	Vitamin K2 (75µg...	16001559	40	from supplier
Hepart	Vitamin K2 (75µg...	18001656	60	from supplier
Hepart	Vitamin K2-Ca-Mn...	16000403	60	from supplier
Hepart	Vitamin K2-Ca-Mn...	15000836	60	from supplier
Hepart	Vitamin K2-Ca-Mn...	15000836	60	from supplier
Hepart	Walnut leaf extr...	15001009	40	from supplier
Hepart	Walnut leaf extr...	19000343/0	40	from supplier
Hepart	Walnut leaf extr...	16001149	40	from supplier
Hepart	Walnut leaf extr...	16001149	40	from supplier
Hepart	Wild berry (H1417)	18000975/0	40	from supplier
Hepart	Wild berry (H1417)	18001313	40	from supplier
Hepart	Yeast organic (H...	17001784/0	40	from supplier
Hepart	Yeast organic (H...	1900527	40	from supplier
Hepart	Yeast organic (H...	16001157	40	from supplier
Hepart	Zinc (H1110)	18000559	40	from supplier
Hepart	Zinc (H1110)	14001803	60	from supplier
Hepart	Zinc (H1110)	16000196	60	from supplier
Hepart	Zinc (H1110)	14001753	60	from supplier
Hepart	Zinc (H1110)	18001012	40	from supplier
Hepart	Zinc (H1110)	16000519	40	from supplier
Hepart	Zinc (H1110)	16000518	40	from supplier
Inresa	Amifampridine	33-16-I-77P3018	40	not necessary
Inresa	Baclofen	206031617/09062017	40	AR-17-FG-018577-01
Inresa	Ketoprofen	F1300183	60	AR-15-FG-009572-01
Inresa	Ketoprofen	F1300182	60	not necessary
Inresa	Maltodextrin	R1503902	60	AR-15-FG-008821-01
Inresa	Maltodextrin	R1600014	60	not necessary
Inresa	Metoprolol tartr...	405005RIIC	40	not necessary
Inresa	Sotalol hydrochl...	369-MC	40	not necessary
Katwijk Chemi...	Ethosuximide	451644	40	not necessary
Klenk	Agarikon	1061A120312	40	pending
Klenk	Angelica root	6085B120614	60	1509356

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Klenk	Angelica root	6085B120614	40	not necessary
Klenk	Aniseed	2700A160511	60	1708049
Klenk	Aniseed	2700A140129	60	not necessary
Klenk	Ash bark, cut	0871B110127	60	not necessary
Klenk	Asparagus root, ...	6151B110304	40	1403480
Klenk	Birch bark, cut	0561A121121	40	not necessary
Klenk	Bitter fennel	2860A140613	60	AR-15-FG-008160-01
Klenk	Bitter fennel	2880N120418	40	AR-15-FG-002805-02
Klenk	Bitter fennel	2880N120418	40	1412547
Klenk	Black mustard seed	7430A120813	40	AR-15-FG-001050-01
Klenk	Black mustard se...	7434A120712	40	AR-15-FG-006711-01
Klenk	Boldo leaf, cut	1911A100714	40	not necessary
Klenk	Buckwheat seed	7510A111011	60	pending
Klenk	Buckwheat seed	7510A111011	40	from supplier
Klenk	Calendula flower	1241B130321	40	not necessary
Klenk	Cardamom, excort...	2765A130524	60	AR-15-FG-015285-01
Klenk	Celery fruit	7140A131122	60	not necessary
Klenk	Celery fruit	7140A130625	40	1405177
Klenk	Celery root, cut	6091A110110	40	1410177
Klenk	Chelidonium maju...	3731K120904	40	not necessary
Klenk	Chervil, rubbed	3721A100329	40	not necessary
Klenk	Cleavers goosegr...	4731B110809	40	not necessary
Klenk	Clove	1260A150916	40	1708050
Klenk	Clove, powdered ...	1264B141007	60	AR-15-FG-014858-01
Klenk	Comfrey root	6244B130624	40	not necessary
Klenk	Cottonwood bud, ...	3266B150624	60	not necessary
Klenk	Cottonwood bud, ...	3266B130826	40	1410178
Klenk	Cumin, powdered	2824B121001	40	not necessary
Klenk	Dandelion root, ...	6621B151029	60	AR-17-FG-007535-01
Klenk	Dill fruit	2680A170222	60	not necessary
Klenk	Elder berries	3010A120202	40	pending
Klenk	Elder leaf, cut	2321B130221	60	not necessary
Klenk	Elecampane rhizo...	6341B120611	40	not necessary
Klenk	Euphrasy herb, cut	3871B100713	40	not necessary
Klenk	Fenugreek, powde...	7264B130205	40	1408080
Klenk	Fenugreek, powde...	7264A151001	60	AR-16-FG-014493-01
Klenk	Galeopsis herb, ...	3911A090219	40	not necessary
Klenk	Ginger, peeled, ...	6954B140724	60	AR-15-FG-014915-01
Klenk	Ginger, peeled, ...	6954B150325	60	1606616
Klenk	Ginger, peeled, ...	6954B160429	60	1704440
Klenk	Ginger, peeled, ...	6954B170127	60	1708675
Klenk	Goat's rue, cut	3901C110905	40	not necessary
Klenk	Greek mountain h...	4511K130314	40	not necessary
Klenk	Green tea	5100A100611	40	not necessary
Klenk	Java tea, cut	2221B110718	40	not necessary
Klenk	Juniper, powdered	2894B130320	40	1411362
Klenk	Juniper, powdered	2894B150922	60	1704439
Klenk	Knotgrass, cut	4301A120912	40	not necessary
Klenk	Ladys bedstraw ...	3921B101012	40	not necessary
Klenk	Laurel leaf, pow...	2124A110329	60	not necessary
Klenk	Laurel leaf, pow...	2124A110329	40	1411032
Klenk	Maize style, cut	5581A111206	40	not necessary
Klenk	Meadowsweet herb...	4541A120313	60	not necessary
Klenk	Melissa leaf, cut	2161A090630	130	not necessary
Klenk	Melissa leaf, cut	2161A121029	55	not necessary
Klenk	Motherwort herb,...	4071A121211	40	not necessary
Klenk	Myrrh	3355B120628	40	AR-15-FG-000945-01
Klenk	Nettle seed	7470A111117	40	pending
Klenk	Nutmeg	7320A121002	40	not necessary
Klenk	Nutmeg, powdered...	7325A160405	60	9888

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
Klenk	Parsley herb, ru...	4262A091116	60	not necessary
Klenk	Parsley seed	7370A130909	40	1407600
Klenk	Passionflower he...	4231A120329	40	not necessary
Klenk	Peony root, cut	6421B120327	60	not necessary
Klenk	Restharrow root,...	6411B130404	40	not necessary
Klenk	Rose hips seed	7240A120402	60	not necessary
Klenk	Rosemary leaf, cut	2261A131118	60	AR-15-FG-008499-01
Klenk	Ruscus rhizome, ...	7081B110301	60	not necessary
Klenk	Saffron (spice)	0384A121022	40	not necessary
Klenk	Senna leaf, powd...	2334B130606	60	AR-15-FG-015212-01
Klenk	Strawberry leaf,...	2031B110224	40	not necessary
Klenk	Strawberry leaf,...	2031B120925	40	not necessary
Klenk	Sweet fennel	2870A130304	40	1405154
Klenk	Sweet fennel	2870A110412	40	pending
Klenk	Theriac, powdered	13816	60	1705088
Klenk	Tormentil	6914B120112	40	not necessary
Klenk	Verbena herb, cut	4651B120611	40	not necessary
Klenk	Willowherb, smal...	4721A140203	60	not necessary
Kräuter Schulte	Twotoothed achyr...	45053	40	not necessary
Kreatin.de GmbH	Creatine	300331	50	1405615
Lycored corp.	Caromix	CMX-LLCV-011401	40	not necessary
Merck KGaA	Benzoyl peroxide...	S7012241509	60	20130317
Otto Fischar	Talcum	90841/1	40	from supplier
Pharmorgana	L-carnitine tart...	1000160672	40	from supplier
Pharmorgana	L-carnitine tart...	1000170113	40	from supplier
Pharmorgana	Sodium dichloroa...	1000150016	40	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708421	40	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708125	40	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708823	60	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708823	60	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708622	60	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708421	60	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708421	75	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708821	60	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708821	60	from supplier
PhytoComm	(Ku) Xing Ren	G029H0708322	40	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532422	40	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532422	75	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532323	40	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532521	60	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532323	75	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532621	60	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532323	40	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532121	40	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532822	60	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532621	60	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532621	60	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532822	60	from supplier
PhytoComm	(Sheng) Di Huang	G211H0532621	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544622	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544421	40	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544622	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544421	40	from supplier
PhytoComm	(Shi) Chang Pu	G006HS2700N1	60	from supplier
PhytoComm	(Shi) Chang Pu	G006HS2700N1	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544821	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544322	40	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544821	60	from supplier
PhytoComm	(Shi) Chang Pu	G006HS2700N1	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544321	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	(Shi) Chang Pu	G006HS270PH1	60	from supplier
PhytoComm	(Shi) Chang Pu	G006H0544421	75	from supplier
PhytoComm	(Shi) Chang Pu	G006HS270PH1	60	from supplier
PhytoComm	Ba Ji Tian	G165H0445821	60	from supplier
PhytoComm	Ba Ji Tian	G165H0445222	40	from supplier
PhytoComm	Ba Ji Tian	G165H0445421	75	from supplier
PhytoComm	Ba Ji Tian	G165H0445621	60	from supplier
PhytoComm	Ba Ji Tian	G165H0445821	60	from supplier
PhytoComm	Ba Ji Tian	G165H0445621	60	from supplier
PhytoComm	Ba Ji Tian	G165H0445621	60	from supplier
PhytoComm	Ba Ji Tian	G165H0445421	40	from supplier
PhytoComm	Ba Ji Tian	G165H0445621	60	from supplier
PhytoComm	Bai Bian Dou	G093H0510522	60	from supplier
PhytoComm	Bai Bian Dou	G093H0510522	60	from supplier
PhytoComm	Bai Bian Dou	G093HS177QM1	60	from supplier
PhytoComm	Bai Bian Dou	G093HS177QM1	60	from supplier
PhytoComm	Bai Bian Dou	G093H0510221	40	from supplier
PhytoComm	Bai Bian Dou	G093H0510821	60	from supplier
PhytoComm	Bai Bian Dou	G093H0510821	60	from supplier
PhytoComm	Bai Bian Dou	G093H0510221	75	from supplier
PhytoComm	Bai Bian Dou	G093H0510221	40	from supplier
PhytoComm	Bai Bian Dou	G093H0510521	60	from supplier
PhytoComm	Bai He	G139H0607321	75	from supplier
PhytoComm	Bai He	G139H0607421	60	from supplier
PhytoComm	Bai He	G139H0607321	40	from supplier
PhytoComm	Bai He	G139H0607321	40	from supplier
PhytoComm	Bai He	G139H0607622	60	from supplier
PhytoComm	Bai He	G139H0607622	60	from supplier
PhytoComm	Bai He	G139H0607221	20	from supplier
PhytoComm	Bai He	G139H0607421	60	from supplier
PhytoComm	Bai Hua She She ...	G173HS078RW1	60	from supplier
PhytoComm	Bai Hua She She ...	G173H1153522	60	from supplier
PhytoComm	Bai Hua She She ...	G173H1153421	60	from supplier
PhytoComm	Bai Hua She She ...	G173H1153522	60	from supplier
PhytoComm	Bai Hua She She ...	G173H1153421	40	from supplier
PhytoComm	Bai Hua She She ...	G173H1153421	59	from supplier
PhytoComm	Bai Hua She She ...	G173H1153222	40	from supplier
PhytoComm	Bai Hua She She ...	G173H1153421	60	from supplier
PhytoComm	Bai Hua She She ...	G173HS078RW1	60	from supplier
PhytoComm	Bai Mao Gen	G125H0512221	40	from supplier
PhytoComm	Bai Mao Gen	G125H0512421	60	from supplier
PhytoComm	Bai Mao Gen	G125H0512521	60	from supplier
PhytoComm	Bai Mao Gen	G125H0512821	60	from supplier
PhytoComm	Bai Mao Gen	G125H0512421	60	from supplier
PhytoComm	Bai Mao Gen	G125H0512821	60	from supplier
PhytoComm	Bai Mao Gen	G125H0512221	40	from supplier
PhytoComm	Bai Shao Yao	G179HS065RQ1	60	from supplier
PhytoComm	Bai Shao Yao	G179HS065RQ1	59	from supplier
PhytoComm	Bai Shao Yao	G179HS065RW1	60	from supplier
PhytoComm	Bai Shao Yao	G179HS065RW1	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502322	40	from supplier
PhytoComm	Bai Shao Yao	G179H0502323	40	from supplier
PhytoComm	Bai Shao Yao	G179H0502422	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502721	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502621	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502721	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502521	60	from supplier
PhytoComm	Bai Shao Yao	G179HS065RT1	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502522	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502422	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Bai Shao Yao	G179HS065RT1	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502522	60	from supplier
PhytoComm	Bai Shao Yao	G179H0502621	60	from supplier
PhytoComm	Bai Xian Pi	G090HS335RP1	60	from supplier
PhytoComm	Bai Xian Pi	G090HS3355PM1	60	from supplier
PhytoComm	Bai Xian Pi	G090HS3355PM1	60	from supplier
PhytoComm	Bai Xian Pi	G090H0515121	40	from supplier
PhytoComm	Bai Xian Pi	G090HS335RP1	60	from supplier
PhytoComm	Bai Xian Pi	G090H0515121	75	from supplier
PhytoComm	Bai Xian Pi	G090HS335RN1	60	from supplier
PhytoComm	Bai Xian Pi	G090H0515521	60	from supplier
PhytoComm	Bai Xian Pi	G090HS335RN1	60	from supplier
PhytoComm	Bai Zhu	G041H0501721	60	from supplier
PhytoComm	Bai Zhu	G041H0501721	60	from supplier
PhytoComm	Bai Zhu	G041H0501422	40	from supplier
PhytoComm	Bai Zhu	G041H0501321	40	from supplier
PhytoComm	Bai Zhu	G041H0501221	40	from supplier
PhytoComm	Bai Zhu	G041H0501321	40	from supplier
PhytoComm	Bai Zhu	G041H0501321	75	from supplier
PhytoComm	Bai Zhu	G041HS064SG1	60	from supplier
PhytoComm	Bai Zhu	G041H0501422	75	from supplier
PhytoComm	Bai Zhu	G041H0501822	60	from supplier
PhytoComm	Bai Zhu	G041H0501822	60	from supplier
PhytoComm	Bai Zhu	G041HS064SG1	60	from supplier
PhytoComm	Bai Zhu	G041H0501621	60	from supplier
PhytoComm	Bai Zi Ren	G052F110315	40	from supplier
PhytoComm	Bai Zi Ren	G052HS187SK1	60	from supplier
PhytoComm	Bai Zi Ren	G052HS187SK1	60	from supplier
PhytoComm	Bai Zi Ren	G052HS187	60	from supplier
PhytoComm	Bai Zi Ren	G052HS187PG1	60	from supplier
PhytoComm	Bai Zi Ren	G052F110320	75	from supplier
PhytoComm	Bai Zi Ren	G052F110320	40	from supplier
PhytoComm	Bai Zi Ren	G052F110320	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526521	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526521	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191HS090SH1	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526322	40	from supplier
PhytoComm	Ban Xia (Jiang)	G191HS090QL1	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526322	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191HS090QL1	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526322	40	from supplier
PhytoComm	Ban Xia (Jiang)	G191HS090SH1	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526421	60	from supplier
PhytoComm	Ban Xia (Jiang)	G191H0526322	75	from supplier
PhytoComm	Ban Zhi Lian	G224H0528522	60	from supplier
PhytoComm	Ban Zhi Lian	G224H0528522	60	from supplier
PhytoComm	Ban Zhi Lian	G224H0528221	40	from supplier
PhytoComm	Ban Zhi Lian	G224H0528521	60	from supplier
PhytoComm	Bu Gu Zhi	G205H1247621	60	from supplier
PhytoComm	Bu Gu Zhi	G205H1247621	60	from supplier
PhytoComm	Bu Gu Zhi	G205H1247421	75	from supplier
PhytoComm	Bu Gu Zhi	G205H1247822	60	from supplier
PhytoComm	Bu Gu Zhi	G205H1247822	60	from supplier
PhytoComm	Bu Gu Zhi	G205H1247222	40	from supplier
PhytoComm	Bu Gu Zhi	G205H1247421	60	from supplier
PhytoComm	Bu Gu Zhi	G205H1247421	40	from supplier
PhytoComm	Cang Er Zi	G251H1401521	60	from supplier
PhytoComm	Cang Er Zi	G251H1401622	60	from supplier
PhytoComm	Cang Er Zi	G251H1401622	60	from supplier
PhytoComm	Cang Er Zi	G251H1401322	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Cang Er Zi	G251H1401622	60	from supplier
PhytoComm	Cang Er Zi	H1401821	30	from supplier
PhytoComm	Cang Er Zi	G251H1401421	75	from supplier
PhytoComm	Cang Er Zi	G251H1401622	60	from supplier
PhytoComm	Cang Er Zi	G251H1401421	40	from supplier
PhytoComm	Cang Er Zi	G251H1401421	40	from supplier
PhytoComm	Cang Er Zi	G251H1401821	60	from supplier
PhytoComm	Cang Er Zi	G251H1401821	60	from supplier
PhytoComm	Cang Zhu	G042H1402322	75	from supplier
PhytoComm	Cang Zhu	G042H1402221	40	from supplier
PhytoComm	Cang Zhu	G042H1402321	40	from supplier
PhytoComm	Cang Zhu	G042H1402621	60	from supplier
PhytoComm	Cang Zhu	G042H1402621	60	from supplier
PhytoComm	Cang Zhu	G042H1402322	40	from supplier
PhytoComm	Cang Zhu	G042HS301SG1	60	from supplier
PhytoComm	Cang Zhu	G042H1402821	60	from supplier
PhytoComm	Cang Zhu	G042H1402521	60	from supplier
PhytoComm	Cang Zhu	G042H1402421	60	from supplier
PhytoComm	Cang Zhu	G042H1402521	60	from supplier
PhytoComm	Cang Zhu	G042H1402322	40	from supplier
PhytoComm	Cang Zhu	G042H1402821	60	from supplier
PhytoComm	Cang Zhu	G042H1402421	60	from supplier
PhytoComm	Cang Zhu	G042HS301SG1	60	from supplier
PhytoComm	Chai Hu	G055H1055721	60	from supplier
PhytoComm	Chai Hu	G055HS233SK1	60	from supplier
PhytoComm	Chai Hu	G055HS233SK1	60	from supplier
PhytoComm	Chai Hu	G055H1055721	60	from supplier
PhytoComm	Chai Hu	G055H1055623	60	from supplier
PhytoComm	Chai Hu	G055H1055623	60	from supplier
PhytoComm	Chai Hu	G055H1055821	60	from supplier
PhytoComm	Chai Hu	G055H1055422	75	from supplier
PhytoComm	Chai Hu	G055H1055322	40	from supplier
PhytoComm	Chai Hu	G055H1055821	60	from supplier
PhytoComm	Chai Hu	G055H1055422	40	from supplier
PhytoComm	Che Qian Zi	G192H0723421	40	from supplier
PhytoComm	Che Qian Zi	G192HS120SH1	60	from supplier
PhytoComm	Che Qian Zi	G192HS120SH1	60	from supplier
PhytoComm	Che Qian Zi	G192H0723323	75	from supplier
PhytoComm	Che Qian Zi	G192H0723321	40	from supplier
PhytoComm	Che Qian Zi	G192H0723621	60	from supplier
PhytoComm	Che Qian Zi	G192H0723421	60	from supplier
PhytoComm	Che Qian Zi	G192H0723421	75	from supplier
PhytoComm	Che Qian Zi	G192H0723322	40	from supplier
PhytoComm	Che Qian Zi	G192H0723323	40	from supplier
PhytoComm	Che Qian Zi	G192H0723621	60	from supplier
PhytoComm	Chen Pi	G073HS254SG1	60	from supplier
PhytoComm	Chen Pi	G073HS254RR1	60	from supplier
PhytoComm	Chen Pi	G073HS254RR1	60	from supplier
PhytoComm	Chen Pi	G073H1128421	40	from supplier
PhytoComm	Chen Pi	G073H1128421	60	from supplier
PhytoComm	Chen Pi	G073H1128222	40	from supplier
PhytoComm	Chen Pi	G073H1128421	75	from supplier
PhytoComm	Chen Pi	G073H1128621	60	from supplier
PhytoComm	Chen Pi	G073H1128221	20	from supplier
PhytoComm	Chen Pi	G073H1128821	60	from supplier
PhytoComm	Chen Pi	G073H1128821	60	from supplier
PhytoComm	Chen Pi	G073HS254SG1	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705422	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705621	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705321	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Chi Shao (Yao)	G180H0705322	40	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705522	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705621	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705322	75	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705522	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705322	40	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705821	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705821	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705822	60	from supplier
PhytoComm	Chi Shao (Yao)	G180H0705822	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306821	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306821	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306621	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306621	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306621	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306621	60	from supplier
PhytoComm	Chuan Lian Zi	G157H0306221	40	from supplier
PhytoComm	Chuan Lian Zi	G157H0306421	60	from supplier
PhytoComm	Chuan Xiong	G136HS031SH1	60	from supplier
PhytoComm	Chuan Xiong	G136HS031SH1	60	from supplier
PhytoComm	Chuan Xiong	G136H0305421	40	from supplier
PhytoComm	Chuan Xiong	G136H0305321	40	from supplier
PhytoComm	Chuan Xiong	G136H0305421	75	from supplier
PhytoComm	Chuan Xiong	G136H0305321	75	from supplier
PhytoComm	Chuan Xiong	G136HS031RR2	60	from supplier
PhytoComm	Chuan Xiong	G136H0305621	60	from supplier
PhytoComm	Chuan Xiong	G136HS031RR2	60	from supplier
PhytoComm	Chuan Xiong	G136H0305621	60	from supplier
PhytoComm	Chuan Xiong	G136H0305622	60	from supplier
PhytoComm	Chuan Xiong	G136H0305622	60	from supplier
PhytoComm	Chuang Mu Xiang	G219H0309322	75	from supplier
PhytoComm	Chuang Mu Xiang	G219H0309521	60	from supplier
PhytoComm	Chuang Mu Xiang	G219HS046SK1	60	from supplier
PhytoComm	Chuang Mu Xiang	G219HS046SK1	60	from supplier
PhytoComm	Chuang Mu Xiang	G219H0309321	40	from supplier
PhytoComm	Chuang Mu Xiang	G219H0309322	40	from supplier
PhytoComm	Chuang Mu Xiang	G219H0309322	40	from supplier
PhytoComm	Ci Wu Jia	G302H0856521	60	from supplier
PhytoComm	Ci Wu Jia	G302H0856322	75	from supplier
PhytoComm	Ci Wu Jia	G302H0856322	40	from supplier
PhytoComm	Ci Wu Jia	G302H0856521	60	from supplier
PhytoComm	Ci Wu Jia	G302H0856322	40	from supplier
PhytoComm	Ci Wu Jia	G302HS422SH1	60	from supplier
PhytoComm	Ci Wu Jia	G302HS422SH1	60	from supplier
PhytoComm	Ci Wu Jia	G302H0856621	60	from supplier
PhytoComm	Ci Wu Jia	G302H0856621	60	from supplier
PhytoComm	Ci Wu Jia	G302H0856421	60	from supplier
PhytoComm	Da Huang	G212H0320622	60	from supplier
PhytoComm	Da Huang	G212H0320622	60	from supplier
PhytoComm	Da Huang	G212H0320622	60	from supplier
PhytoComm	Da Huang	G212H0320521	60	from supplier
PhytoComm	Da Huang	G212H0320322	40	from supplier
PhytoComm	Da Huang	G212H0320521	60	from supplier
PhytoComm	Da Huang	G212H0320822	60	from supplier
PhytoComm	Da Huang	G212H0320322	40	from supplier
PhytoComm	Da Huang	G212H0320622	60	from supplier
PhytoComm	Da Huang	G212H0320521	60	from supplier
PhytoComm	Da Huang	G212H0320321	41	from supplier
PhytoComm	Da Huang	G212H0320521	60	from supplier
PhytoComm	Da Huang	G212H0320822	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Da Huang	G212H0320322	75	from supplier
PhytoComm	Da Zao	G129H0316621	60	from supplier
PhytoComm	Da Zao	G129H0316321	60	from supplier
PhytoComm	Da Zao	G129H0316521	60	from supplier
PhytoComm	Da Zao	G129H0316321	40	from supplier
PhytoComm	Da Zao	G129H0316321	75	from supplier
PhytoComm	Da Zao	G129H0316821	60	from supplier
PhytoComm	Da Zao	G129H0316821	60	from supplier
PhytoComm	Da Zao	G129HS019RT1	60	from supplier
PhytoComm	Da Zao	G129H0316621	60	from supplier
PhytoComm	Da Zao	G129HS019RT1	60	from supplier
PhytoComm	Dan Dou Chi	G232HS239SK1	60	from supplier
PhytoComm	Dan Dou Chi	G232HS239	60	from supplier
PhytoComm	Dan Dou Chi	G232HS239SK1	60	from supplier
PhytoComm	Dan Dou Chi	G232H1106321	40	from supplier
PhytoComm	Dan Dou Chi	G232H1106621	55	from supplier
PhytoComm	Dan Dou Chi	G232H1106621	60	from supplier
PhytoComm	Dan Dou Chi	G232H1106321	75	from supplier
PhytoComm	Dan Dou Chi	G232H1106321	40	from supplier
PhytoComm	Dan Dou Chi	G232H1106621	60	from supplier
PhytoComm	Dan Dou Chi	G232H1106621	60	from supplier
PhytoComm	Dan Shen	G214H0441322	40	from supplier
PhytoComm	Dan Shen	G214H0441324	40	from supplier
PhytoComm	Dan Shen	G214H0441324	75	from supplier
PhytoComm	Dan Shen	G214H0441821	60	from supplier
PhytoComm	Dan Shen	G214HS059SG1	60	from supplier
PhytoComm	Dan Shen	G214H0441522	60	from supplier
PhytoComm	Dan Shen	G214H0441522	60	from supplier
PhytoComm	Dan Shen	G214H0441821	60	from supplier
PhytoComm	Dan Shen	G214H0441422	60	from supplier
PhytoComm	Dan Shen	G214H0441522	60	from supplier
PhytoComm	Dan Shen	G214HS059SG1	60	from supplier
PhytoComm	Dan Shen	G214H0441324	20	from supplier
PhytoComm	Dan Zhu Ye	G046H1105321	75	from supplier
PhytoComm	Dan Zhu Ye	G046H1105321	40	from supplier
PhytoComm	Dan Zhu Ye	G046H1105321	40	from supplier
PhytoComm	Dan Zhu Ye	G046H1105821	60	from supplier
PhytoComm	Dan Zhu Ye	G046H1105821	60	from supplier
PhytoComm	Dan Zhu Ye	G046H1105621	60	from supplier
PhytoComm	Dan Zhu Ye	G046H1105621	60	from supplier
PhytoComm	Dang Gui	G022H1306422	75	from supplier
PhytoComm	Dang Gui	G022H1306121	40	from supplier
PhytoComm	Dang Gui	G022HS290SK1	60	from supplier
PhytoComm	Dang Gui	G022HS290SK1	60	from supplier
PhytoComm	Dang Gui	G022H1306422	40	from supplier
PhytoComm	Dang Gui	G022H1306721	60	from supplier
PhytoComm	Dang Gui	G022H1306621	60	from supplier
PhytoComm	Dang Gui	G022H1306321	40	from supplier
PhytoComm	Dang Gui	G022H1306721	60	from supplier
PhytoComm	Dang Gui	G022H1306023	40	from supplier
PhytoComm	Dang Gui	G022H1306822	60	from supplier
PhytoComm	Dang Gui	G022H1306822	60	from supplier
PhytoComm	Dang Gui Wei	G312HS291SK1	60	from supplier
PhytoComm	Dang Gui Wei	G312H1307621	60	from supplier
PhytoComm	Dang Gui Wei	G312H1307522	60	from supplier
PhytoComm	Dang Gui Wei	G312H1307621	60	from supplier
PhytoComm	Dang Gui Wei	G312H1307221	60	from supplier
PhytoComm	Dang Gui Wei	G312HS291SK1	60	from supplier
PhytoComm	Dang Shen	G077H2001621	60	from supplier
PhytoComm	Dang Shen	G077H2001621	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Dang Shen	G077H2001321	40	from supplier
PhytoComm	Dang Shen	G077H2001521	60	from supplier
PhytoComm	Dang Shen	G077H2001821	60	from supplier
PhytoComm	Dang Shen	G077H2001821	60	from supplier
PhytoComm	Dang Shen	G077HS352SG1	60	from supplier
PhytoComm	Dang Shen	G077HS352RQ1	60	from supplier
PhytoComm	Dang Shen	G077HS352RQ1	60	from supplier
PhytoComm	Dang Shen	G077H2001421	75	from supplier
PhytoComm	Dang Shen	G077H2001421	40	from supplier
PhytoComm	Dang Shen	G077HS352SG1	60	from supplier
PhytoComm	Dang Shen	G077H2001222	20	from supplier
PhytoComm	Di Gu Pi	G149H0601421	75	from supplier
PhytoComm	Di Gu Pi	G149H0601421	40	from supplier
PhytoComm	Di Gu Pi	G149H0601622	60	from supplier
PhytoComm	Di Gu Pi	G149HS101SK1	60	from supplier
PhytoComm	Di Gu Pi	G149H0601821	60	from supplier
PhytoComm	Di Gu Pi	G149HS101SK1	60	from supplier
PhytoComm	Di Gu Pi	G149H0601821	60	from supplier
PhytoComm	Di Gu Pi	G149H0601622	60	from supplier
PhytoComm	Dong Gua Zi	G049H0551421	60	from supplier
PhytoComm	Dong Gua Zi	G049H0551421	75	from supplier
PhytoComm	Dong Gua Zi	G049H0551421	40	from supplier
PhytoComm	Dong Gua Zi	G049H0551622	60	from supplier
PhytoComm	Dong Gua Zi	G049H0551622	60	from supplier
PhytoComm	Dong Gua Zi	G049H0551822	60	from supplier
PhytoComm	Dong Gua Zi	G049H0551822	60	from supplier
PhytoComm	Dong Kui Zi	G001H0834321	38	from supplier
PhytoComm	Dong Kui Zi	G001H0834321	60	from supplier
PhytoComm	Dong Kui Zi	G001H0834821	60	from supplier
PhytoComm	Dong Kui Zi	G001H0834821	60	from supplier
PhytoComm	Dong Kui Zi	G001H0834721	60	from supplier
PhytoComm	Dong Kui Zi	G001H0834721	60	from supplier
PhytoComm	Du Huo	G021H1610221	40	from supplier
PhytoComm	Du Huo	G021H1610421	75	from supplier
PhytoComm	Du Huo	G021H1610821	60	from supplier
PhytoComm	Du Huo	G021H1610821	60	from supplier
PhytoComm	Du Huo	G021H1610421	40	from supplier
PhytoComm	Du Zhong	G101HS133SH1	60	from supplier
PhytoComm	Du Zhong	G101HS133SH1	60	from supplier
PhytoComm	Du Zhong	G101H0736421	75	from supplier
PhytoComm	Du Zhong	G101H0736221	40	from supplier
PhytoComm	Du Zhong	G101H0736522	60	from supplier
PhytoComm	Du Zhong	G101H0736621	60	from supplier
PhytoComm	Du Zhong	G101H0736821	60	from supplier
PhytoComm	Du Zhong	G101H0736622	60	from supplier
PhytoComm	Du Zhong	G101H0736821	60	from supplier
PhytoComm	Du Zhong	G101H0736421	40	from supplier
PhytoComm	Du Zhong	G101H0736321	40	from supplier
PhytoComm	Du Zhong	G101H0736622	60	from supplier
PhytoComm	Du Zhong	G101H0736621	60	from supplier
PhytoComm	E Zhu	G083H1110521	60	from supplier
PhytoComm	E Zhu	G083H1110222	40	from supplier
PhytoComm	E Zhu	G083H1110821	60	from supplier
PhytoComm	E Zhu	G083H1110222	75	from supplier
PhytoComm	E Zhu	G083H1110821	60	from supplier
PhytoComm	E Zhu	G083H1110222	40	from supplier
PhytoComm	Fang Feng	G217H0728521	60	from supplier
PhytoComm	Fang Feng	G217H0728421	40	from supplier
PhytoComm	Fang Feng	G217H0728521	60	from supplier
PhytoComm	Fang Feng	G217H0728821	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Fang Feng	G217H0728621	60	from supplier
PhytoComm	Fang Feng	G217H0728621	60	from supplier
PhytoComm	Fang Feng	G217H0728521	60	from supplier
PhytoComm	Fang Feng	G217H0728821	60	from supplier
PhytoComm	Fang Feng	G217H0728421	75	from supplier
PhytoComm	Fang Ji	G236HS122QM1	60	from supplier
PhytoComm	Fang Ji	G236H0727321	40	from supplier
PhytoComm	Fang Ji	G236H0727321	40	from supplier
PhytoComm	Fang Ji	G236H0727123	40	from supplier
PhytoComm	Fang Ji	G236HS122QM1	60	from supplier
PhytoComm	Fang Ji	G236H0727321	75	from supplier
PhytoComm	Fo Shou	G307HS159RQ1	60	from supplier
PhytoComm	Fo Shou	G307H0888321	30	from supplier
PhytoComm	Fo Shou	G307H0888321	75	from supplier
PhytoComm	Fo Shou	G307HS159RQ1	60	from supplier
PhytoComm	Fo Shou	G307H0888321	40	from supplier
PhytoComm	Fo Shou	G307HS159PM1	60	from supplier
PhytoComm	Fo Shou	G307HS159PM1	60	from supplier
PhytoComm	Fu Ling	G200H1019623	60	from supplier
PhytoComm	Fu Ling	G200HS211SH1	60	from supplier
PhytoComm	Fu Ling	G200H1019421	75	from supplier
PhytoComm	Fu Ling	G200H1019623	60	from supplier
PhytoComm	Fu Ling	G200H1019322	75	from supplier
PhytoComm	Fu Ling	G200H1019622	60	from supplier
PhytoComm	Fu Ling	G200H1019622	60	from supplier
PhytoComm	Fu Ling	G200HS211	60	from supplier
PhytoComm	Fu Ling	G200HS211	60	from supplier
PhytoComm	Fu Ling	G200H1019321	40	from supplier
PhytoComm	Fu Ling	G200H1019421	40	from supplier
PhytoComm	Fu Ling	G200H1019322	40	from supplier
PhytoComm	Fu Ling	G200HS211RT1	60	from supplier
PhytoComm	Fu Ling	G200H1019522	60	from supplier
PhytoComm	Fu Ling	G200HS211RT1	60	from supplier
PhytoComm	Fu Ling	G200H1019322	40	from supplier
PhytoComm	Fu Ling	G200HS211SH1	60	from supplier
PhytoComm	Fu Shen	G201H1021422	40	from supplier
PhytoComm	Fu Shen	G201HS212	60	from supplier
PhytoComm	Fu Shen	G201H1021321	41	from supplier
PhytoComm	Fu Shen	G201HS212SH1	60	from supplier
PhytoComm	Fu Shen	G201HS212SH1	60	from supplier
PhytoComm	Fu Shen	G201H1021422	60	from supplier
PhytoComm	Fu Shen	G201H1021422	52	from supplier
PhytoComm	Fu Xiao Mai	G244HS021SK1	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS021QM1	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS021QM1	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS0210V1	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS021RT1	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS021RT1	60	from supplier
PhytoComm	Fu Xiao Mai	G244H1078421	75	from supplier
PhytoComm	Fu Xiao Mai	G244H1078421	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS0210V1	60	from supplier
PhytoComm	Fu Xiao Mai	G244HS021SK1	60	from supplier
PhytoComm	Fu Xiao Mai	G244H1078421	40	from supplier
PhytoComm	Fu Zi	G005H0801723	60	from supplier
PhytoComm	Fu Zi	G005H0801122	40	from supplier
PhytoComm	Fu Zi	G005H0801723	60	from supplier
PhytoComm	Fu Zi	G005H0801521	60	from supplier
PhytoComm	Fu Zi	G005HS158SP1	60	from supplier
PhytoComm	Fu Zi	G005H0801524	60	from supplier
PhytoComm	Fu Zi	G005H0801123	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Fu Zi	G005H0801123	60	from supplier
PhytoComm	Fu Zi	G005H0801123	75	from supplier
PhytoComm	Fu Zi	G005H0801624	60	from supplier
PhytoComm	Fu Zi	G005H0801624	60	from supplier
PhytoComm	Fu Zi	G005H0801322	40	from supplier
PhytoComm	Fu Zi	G005H0801524	60	from supplier
PhytoComm	Fu Zi	G005HS158SP1	60	from supplier
PhytoComm	Fu Zi	G005H0801524	60	from supplier
PhytoComm	Gan Cao	G119H0521622	60	from supplier
PhytoComm	Gan Cao	G119H0521121	40	from supplier
PhytoComm	Gan Cao	G119H0521422	40	from supplier
PhytoComm	Gan Cao	G119H0521524	60	from supplier
PhytoComm	Gan Cao	G119H0521823	60	from supplier
PhytoComm	Gan Cao	G119H0521823	60	from supplier
PhytoComm	Gan Cao	G119H0521524	60	from supplier
PhytoComm	Gan Cao	G119H0521323	40	from supplier
PhytoComm	Gan Cao	G119H0521121	75	from supplier
PhytoComm	Gan Cao	G119H0521822	60	from supplier
PhytoComm	Gan Cao	G119H0521522	60	from supplier
PhytoComm	Gan Cao	G119H0521422	75	from supplier
PhytoComm	Gan Cao	G119H0521822	60	from supplier
PhytoComm	Gan Cao	G119H0521422	60	from supplier
PhytoComm	Gan Cao	G119H0521622	60	from supplier
PhytoComm	Gan Jiang	G252H1145521	60	from supplier
PhytoComm	Gan Jiang	G252H1145821	60	from supplier
PhytoComm	Gan Jiang	G252H1145421	40	from supplier
PhytoComm	Gan Jiang	G252H1145821	60	from supplier
PhytoComm	Gan Jiang	G252H1145421	75	from supplier
PhytoComm	Gan Jiang	G252H1145221	40	from supplier
PhytoComm	Gan Jiang	G252H1145621	60	from supplier
PhytoComm	Gan Jiang	G252H1145321	40	from supplier
PhytoComm	Gan Jiang	G252H1145622	60	from supplier
PhytoComm	Gan Jiang	G252H1145622	60	from supplier
PhytoComm	Gan Jiang	G252H1145621	60	from supplier
PhytoComm	Ge Gen	G206H1301422	75	from supplier
PhytoComm	Ge Gen	G206H1301221	40	from supplier
PhytoComm	Ge Gen	G206H1301422	60	from supplier
PhytoComm	Ge Gen	G206H1301822	60	from supplier
PhytoComm	Ge Gen	G206H1301822	60	from supplier
PhytoComm	Ge Gen	G206H1301621	60	from supplier
PhytoComm	Ge Gen	G206H1301323	40	from supplier
PhytoComm	Ge Gen	G206H1301422	40	from supplier
PhytoComm	Ge Gen	G206H1301621	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913422	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913622	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913622	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913821	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913221	100	from supplier
PhytoComm	Gou Qi Zi	G150H0913221	40	from supplier
PhytoComm	Gou Qi Zi	G150H0913221	40	from supplier
PhytoComm	Gou Qi Zi	G150HS171SK1	60	from supplier
PhytoComm	Gou Qi Zi	G150HS171SK1	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913821	60	from supplier
PhytoComm	Gou Qi Zi	G150H0913621	60	from supplier
PhytoComm	Gou Teng	G247H1210621	60	from supplier
PhytoComm	Gou Teng	G247H1210621	60	from supplier
PhytoComm	Gou Teng	G247H1210321	40	from supplier
PhytoComm	Gou Teng	G247H1210322	75	from supplier
PhytoComm	Gou Teng	G247HS288SG1	60	from supplier
PhytoComm	Gou Teng	G247HS288SG1	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Gou Teng	G247H1210822	60	from supplier
PhytoComm	Gou Teng	G247H1210822	60	from supplier
PhytoComm	Gou Teng	G247H1210322	40	from supplier
PhytoComm	Gou Teng	G247H1210322	40	from supplier
PhytoComm	Gou Teng	G247H1210422	60	from supplier
PhytoComm	Gua Lou	G241H1096523	60	from supplier
PhytoComm	Gua Lou	G241H1096321	40	from supplier
PhytoComm	Gua Lou	G241H1096621	60	from supplier
PhytoComm	Gua Lou	G241H1096621	60	from supplier
PhytoComm	Gua Lou	G241H1096521	60	from supplier
PhytoComm	Gua Lou	G241H1096821	60	from supplier
PhytoComm	Gua Lou	G241HS202SK1	60	from supplier
PhytoComm	Gua Lou	G241H1096523	60	from supplier
PhytoComm	Gua Lou	G241H1096523	59	from supplier
PhytoComm	Gua Lou	G241H1096821	60	from supplier
PhytoComm	Gua Lou	G241H1096523	60	from supplier
PhytoComm	Gua Lou	G241HS202SK1	60	from supplier
PhytoComm	Gua Lou	G241H1096322	40	from supplier
PhytoComm	Guang Huo Xiang	G007H2005621	60	from supplier
PhytoComm	Guang Huo Xiang	G007H2005423	60	from supplier
PhytoComm	Guang Huo Xiang	G007H2005521	60	from supplier
PhytoComm	Guang Huo Xiang	G007H2005822	60	from supplier
PhytoComm	Guang Huo Xiang	G007H2005423	40	from supplier
PhytoComm	Guang Huo Xiang	G007H2005822	60	from supplier
PhytoComm	Guang Huo Xiang	G007H2005621	60	from supplier
PhytoComm	Guang Huo Xiang	G007H2005423	75	from supplier
PhytoComm	Gui Zhi	G070H1001422	75	from supplier
PhytoComm	Gui Zhi	G070H1001422	60	from supplier
PhytoComm	Gui Zhi	G070H1001624	60	from supplier
PhytoComm	Gui Zhi	G070H1001624	60	from supplier
PhytoComm	Gui Zhi	G070H1001621	60	from supplier
PhytoComm	Gui Zhi	G070H1001621	60	from supplier
PhytoComm	Gui Zhi	G070H1001621	60	from supplier
PhytoComm	Gui Zhi	G070H1001621	60	from supplier
PhytoComm	Gui Zhi	G070H1001322	40	from supplier
PhytoComm	Gui Zhi	G070H1001621	60	from supplier
PhytoComm	Gui Zhi	G070HS205SH1	60	from supplier
PhytoComm	Gui Zhi	G070H1001621	60	from supplier
PhytoComm	Gui Zhi	G070H1001523	60	from supplier
PhytoComm	Gui Zhi	G070H1001422	40	from supplier
PhytoComm	Gui Zhi	G070HS205SH1	60	from supplier
PhytoComm	Gui Zhi	G070H1001221	40	from supplier
PhytoComm	Gui Zhi	G070H1001821	60	from supplier
PhytoComm	Gui Zhi	G070H1001821	60	from supplier
PhytoComm	Gui Zhi	G070H1001822	60	from supplier
PhytoComm	Gui Zhi	G070H1001822	60	from supplier
PhytoComm	Han Lian Cao	G095H0770321	40	from supplier
PhytoComm	Han Lian Cao	G095H0770622	60	from supplier
PhytoComm	Han Lian Cao	G095H0770622	60	from supplier
PhytoComm	Han Lian Cao	G095HS141RW1	60	from supplier
PhytoComm	Han Lian Cao	G095HS141RW1	60	from supplier
PhytoComm	Han Lian Cao	G095H0770521	60	from supplier
PhytoComm	Han Lian Cao	G095H0770421	40	from supplier
PhytoComm	Han Lian Cao	G095H0770321	40	from supplier
PhytoComm	Han Lian Cao	G095H0770421	75	from supplier
PhytoComm	Han Lian Cao	G095H0770321	75	from supplier
PhytoComm	Han Lian Cao	G095H0770421	60	from supplier
PhytoComm	He Huan Pi	G010H0635822	60	from supplier
PhytoComm	He Huan Pi	G010H0635822	60	from supplier
PhytoComm	He Huan Pi	G010H0635523	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	He Huan Pi	G010H0635321	40	from supplier
PhytoComm	He Huan Pi	G010H0635622	60	from supplier
PhytoComm	He Huan Pi	G010H0635622	60	from supplier
PhytoComm	He Huan Pi	G010H0635423	60	from supplier
PhytoComm	He Huan Pi	G010H0635422	40	from supplier
PhytoComm	He Huan Pi	G010HS114SK1	60	from supplier
PhytoComm	He Huan Pi	G010HS114SK1	60	from supplier
PhytoComm	He Huan Pi	G010H0635422	75	from supplier
PhytoComm	He Huan Pi	G010H0635423	60	from supplier
PhytoComm	He Huan Pi	G010H0635523	60	from supplier
PhytoComm	He Shou Wu	G198H0718421	60	from supplier
PhytoComm	He Shou Wu	G198H0718421	60	from supplier
PhytoComm	He Shou Wu	G198H0718421	60	from supplier
PhytoComm	He Shou Wu	G198H0718421	40	from supplier
PhytoComm	He Shou Wu	G198H0718222	40	from supplier
PhytoComm	He Shou Wu	G198H0718421	75	from supplier
PhytoComm	He Shou Wu	G198H0718821	60	from supplier
PhytoComm	He Shou Wu	G198H0718821	60	from supplier
PhytoComm	Hong Hua	G057H0901621	60	from supplier
PhytoComm	Hong Hua	G057H0901322	75	from supplier
PhytoComm	Hong Hua	G057H0901322	40	from supplier
PhytoComm	Hong Hua	G057H0901621	60	from supplier
PhytoComm	Hong Hua	G057H0901322	40	from supplier
PhytoComm	Hong Hua	G057H0901221	20	from supplier
PhytoComm	Hong Hua	G057H0901622	60	from supplier
PhytoComm	Hong Hua	G057H0901622	60	from supplier
PhytoComm	Hong Hua	G057HS173SG1	60	from supplier
PhytoComm	Hong Hua	G057H0901221	40	from supplier
PhytoComm	Hong Hua	G057HS173SG1	60	from supplier
PhytoComm	Hong Jing Tian	G237H0906321	75	from supplier
PhytoComm	Hong Jing Tian	G237H0906321	40	from supplier
PhytoComm	Hong Jing Tian	G237H0906621	60	from supplier
PhytoComm	Hong Jing Tian	G237H0906821	60	from supplier
PhytoComm	Hong Jing Tian	G237H0906421	40	from supplier
PhytoComm	Hong Jing Tian	G237H0906421	59	from supplier
PhytoComm	Hong Jing Tian	G237H0906321	40	from supplier
PhytoComm	Hong Jing Tian	G237H0906621	60	from supplier
PhytoComm	Hong Jing Tian	G237H0906421	75	from supplier
PhytoComm	Hong Jing Tian	G237H0906821	60	from supplier
PhytoComm	Hong Jing Tian	G237HS429SK1	60	from supplier
PhytoComm	Hong Jing Tian	G237HS429SK1	60	from supplier
PhytoComm	Hou Po	G154H0917421	100	from supplier
PhytoComm	Hou Po	G154H0917521	60	from supplier
PhytoComm	Hou Po	G154H0917721	60	from supplier
PhytoComm	Hou Po	G154H0917423	60	from supplier
PhytoComm	Hou Po	G154H0917822	60	from supplier
PhytoComm	Hou Po	G154H0917721	60	from supplier
PhytoComm	Hou Po	G154H0917222	40	from supplier
PhytoComm	Hou Po	G154H0917421	40	from supplier
PhytoComm	Hou Po	G154H0917521	60	from supplier
PhytoComm	Hou Po	G154H0917822	60	from supplier
PhytoComm	Hou Po	G154HS184SH1	60	from supplier
PhytoComm	Hou Po	G154HS184SH1	60	from supplier
PhytoComm	Hu Zhang	G197HS373SK1	60	from supplier
PhytoComm	Hu Zhang	G197H0851621	60	from supplier
PhytoComm	Hu Zhang	G197HS373SK1	60	from supplier
PhytoComm	Hu Zhang	G197H0851321	40	from supplier
PhytoComm	Hu Zhang	G197H0851321	40	from supplier
PhytoComm	Hu Zhang	G197H0851521	60	from supplier
PhytoComm	Hu Zhang	G197H0851321	75	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Hu Zhang	G197H0851621	60	from supplier
PhytoComm	Hua Shi	G238H1332521	60	from supplier
PhytoComm	Hua Shi	G238H1332621	45	from supplier
PhytoComm	Hua Shi	G238HS289RM1	60	from supplier
PhytoComm	Hua Shi	G238H1332621	48	from supplier
PhytoComm	Hua Shi	G238HS289RM1	60	from supplier
PhytoComm	Hua Shi	G238H1332321	40	from supplier
PhytoComm	Hua Shi	G238HS289SK1	60	from supplier
PhytoComm	Hua Shi	G238HS289SK1	60	from supplier
PhytoComm	Huang Bai	G188HS262SH1	60	from supplier
PhytoComm	Huang Bai	G188HS262SH1	60	from supplier
PhytoComm	Huang Bai	G188H1203322	40	from supplier
PhytoComm	Huang Bo	G188H1203622	60	from supplier
PhytoComm	Huang Bo	G188H1203622	60	from supplier
PhytoComm	Huang Bo	G188H1203321	75	from supplier
PhytoComm	Huang Bo	G188H1203321	40	from supplier
PhytoComm	Huang Bo	G188H1203823	60	from supplier
PhytoComm	Huang Bo	G188H1203823	60	from supplier
PhytoComm	Huang Bo	G188H1203321	40	from supplier
PhytoComm	Huang Jing	FA1206822	40	from supplier
PhytoComm	Huang Jing	G196H1206521	60	from supplier
PhytoComm	Huang Jing	G196H1206321	75	from supplier
PhytoComm	Huang Jing	G196H1206721	60	from supplier
PhytoComm	Huang Jing	G196H1206522	60	from supplier
PhytoComm	Huang Jing	G196H1206721	60	from supplier
PhytoComm	Huang Jing	G196H1206321	40	from supplier
PhytoComm	Huang Jing	G196H1206321	40	from supplier
PhytoComm	Huang Jing	G196H1206522	60	from supplier
PhytoComm	Huang Lian	G079HS261SK1	60	from supplier
PhytoComm	Huang Lian	G079HS261RM1	60	from supplier
PhytoComm	Huang Lian	G079HS261RM1	60	from supplier
PhytoComm	Huang Lian	G079H1205621	60	from supplier
PhytoComm	Huang Lian	G079H1205421	75	from supplier
PhytoComm	Huang Lian	G079H1205321	40	from supplier
PhytoComm	Huang Lian	G079H1205421	40	from supplier
PhytoComm	Huang Lian	G079H1205521	60	from supplier
PhytoComm	Huang Lian	G079H1205622	60	from supplier
PhytoComm	Huang Lian	G079HS261SK1	60	from supplier
PhytoComm	Huang Lian	G079H1205622	60	from supplier
PhytoComm	Huang Lian	FA1205702	40	from supplier
PhytoComm	Huang Lian	G079H1205621	60	from supplier
PhytoComm	Huang Lian	G079H1205421	60	from supplier
PhytoComm	Huang Qi	G040H1202321	40	from supplier
PhytoComm	Huang Qi	G040H1202621	60	from supplier
PhytoComm	Huang Qi	G040H1202621	60	from supplier
PhytoComm	Huang Qi	G040H1202222	40	from supplier
PhytoComm	Huang Qi	G040H1202822	60	from supplier
PhytoComm	Huang Qi	G040H1202420	40	from supplier
PhytoComm	Huang Qi	G040H1202423	40	from supplier
PhytoComm	Huang Qi	G040H1202423	75	from supplier
PhytoComm	Huang Qi	G040H1202822	60	from supplier
PhytoComm	Huang Qi	G040H1202621	60	from supplier
PhytoComm	Huang Qi	G040H1202823	60	from supplier
PhytoComm	Huang Qi	G040H1202823	60	from supplier
PhytoComm	Huang Qi	G040H1202323	40	from supplier
PhytoComm	Huang Qi	G040H1202721	60	from supplier
PhytoComm	Huang Qi	G040H1202521	60	from supplier
PhytoComm	Huang Qi	G040H1202621	60	from supplier
PhytoComm	Huang Qi	G040H1202420	75	from supplier
PhytoComm	Huang Qi	G040H1202721	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Huang Qin	G225HS259SM1	60	from supplier
PhytoComm	Huang Qin	G225H1204322	75	from supplier
PhytoComm	Huang Qin	G225H1204422	40	from supplier
PhytoComm	Huang Qin	G225HS259RM1	60	from supplier
PhytoComm	Huang Qin	G225H1204522	60	from supplier
PhytoComm	Huang Qin	G225HS259RM1	60	from supplier
PhytoComm	Huang Qin	G225H1204322	40	from supplier
PhytoComm	Huang Qin	G225H1204622	60	from supplier
PhytoComm	Huang Qin	G225H1204622	60	from supplier
PhytoComm	Huang Qin	G225H1204422	75	from supplier
PhytoComm	Huang Qin	G225H1204422	40	from supplier
PhytoComm	Huang Qin	G225H1204322	35	from supplier
PhytoComm	Huang Qin	G225HS259SM1	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449321	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449321	75	from supplier
PhytoComm	Huo Ma Ren	G141HS061	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449823	60	from supplier
PhytoComm	Huo Ma Ren	G141HS061SP1	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449823	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449621	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449621	60	from supplier
PhytoComm	Huo Ma Ren	G141HS061SP1	60	from supplier
PhytoComm	Huo Ma Ren	G141H0449321	40	from supplier
PhytoComm	Ji Li	G240H1410421	40	from supplier
PhytoComm	Ji Li	G240H1410321	40	from supplier
PhytoComm	Ji Li	G240H1410321	75	from supplier
PhytoComm	Ji Li	G240H1410821	60	from supplier
PhytoComm	Ji Li	G240H1410522	60	from supplier
PhytoComm	Ji Li	G240H1410522	60	from supplier
PhytoComm	Ji Li	G240H1410321	39	from supplier
PhytoComm	Ji Li	G240H1410821	60	from supplier
PhytoComm	Ji Li	G240HS305SK1	60	from supplier
PhytoComm	Ji Li	G240H1410421	75	from supplier
PhytoComm	Ji Li	G240HS305SK1	60	from supplier
PhytoComm	Ji Li	G240H1410221	40	from supplier
PhytoComm	Ji Xue Teng	G159H2101122	40	from supplier
PhytoComm	Ji Xue Teng	G159H2101622	60	from supplier
PhytoComm	Ji Xue Teng	G159H2101622	60	from supplier
PhytoComm	Ji Xue Teng	G159H2101622	60	from supplier
PhytoComm	Ji Xue Teng	G159H2101622	60	from supplier
PhytoComm	Ji Xue Teng	G159HS337SK1	60	from supplier
PhytoComm	Ji Xue Teng	G159HS337SK1	60	from supplier
PhytoComm	Ji Xue Teng	G159H2101423	60	from supplier
PhytoComm	Ji Xue Teng	G159H2101122	75	from supplier
PhytoComm	Jiang Huang	G085H1716422	60	from supplier
PhytoComm	Jiang Huang	G085H1716321	75	from supplier
PhytoComm	Jiang Huang	G085H1716821	60	from supplier
PhytoComm	Jiang Huang	G085H1716721	60	from supplier
PhytoComm	Jiang Huang	G085H1716721	60	from supplier
PhytoComm	Jiang Huang	G085H1716321	40	from supplier
PhytoComm	Jiang Huang	G085H1716321	40	from supplier
PhytoComm	Jiang Huang	G085H1716821	60	from supplier
PhytoComm	Jiao Gu Lan	G065HS424SG1	60	from supplier
PhytoComm	Jiao Gu Lan	G065F-B1276321	40	from supplier
PhytoComm	Jiao Gu Lan	G065FB1276621	60	from supplier
PhytoComm	Jiao Gu Lan	G065FB1276621	60	from supplier
PhytoComm	Jiao Gu Lan	G065FB1276421	40	from supplier
PhytoComm	Jiao Gu Lan	G065HS424RR1	60	from supplier
PhytoComm	Jiao Gu Lan	G065HS424RR1	60	from supplier
PhytoComm	Jiao Gu Lan	G065FB1276421	75	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Jiao Gu Lan	G065HS424SG1	60	from supplier
PhytoComm	Jiao Gu Lan	G065FB1276421	60	from supplier
PhytoComm	Jie Geng	G193H1015621	60	from supplier
PhytoComm	Jie Geng	G193H1015621	60	from supplier
PhytoComm	Jie Geng	G193H1015821	60	from supplier
PhytoComm	Jie Geng	G193H1015322	40	from supplier
PhytoComm	Jie Geng	G193H1015821	60	from supplier
PhytoComm	Jie Geng	G193H1015721	60	from supplier
PhytoComm	Jie Geng	G193H1015721	60	from supplier
PhytoComm	Jie Geng	G193H1015521	60	from supplier
PhytoComm	Jie Geng	G193H1015322	75	from supplier
PhytoComm	Jie Geng	G193H1015221	40	from supplier
PhytoComm	Jie Geng	G193H1015322	60	from supplier
PhytoComm	Jin Qian Cao	G151H0847821	60	from supplier
PhytoComm	Jin Qian Cao	G151HS368	60	from supplier
PhytoComm	Jin Qian Cao	G151H0847421	75	from supplier
PhytoComm	Jin Qian Cao	G151H0847821	60	from supplier
PhytoComm	Jin Qian Cao	G151HS368	60	from supplier
PhytoComm	Jin Qian Cao	G151H0847421	40	from supplier
PhytoComm	Jin Yin Hua	G145H0845621	60	from supplier
PhytoComm	Jin Yin Hua	G145H0845321	40	from supplier
PhytoComm	Jin Yin Hua	G145HS311SK1	60	from supplier
PhytoComm	Jin Yin Hua	G145HS311SK1	60	from supplier
PhytoComm	Jin Yin Hua	G145H0845321	75	from supplier
PhytoComm	Jin Yin Hua	G145H0845621	60	from supplier
PhytoComm	Jin Yin Hua	G145H0845621	60	from supplier
PhytoComm	Jin Yin Hua	G145HS311QG1	60	from supplier
PhytoComm	Jin Yin Hua	G145H0845621	60	from supplier
PhytoComm	Jin Yin Hua	G145H0845321	40	from supplier
PhytoComm	Jin Yin Hua	G145H0845323	60	from supplier
PhytoComm	Jin Yin Hua	G145H0845323	60	from supplier
PhytoComm	Jin Yin Hua	H0845923	30	from supplier
PhytoComm	Jin Yin Hua	G145HS311QG1	60	from supplier
PhytoComm	Jing Jie	G221H1038421	75	from supplier
PhytoComm	Jing Jie	G221HS218	60	from supplier
PhytoComm	Jing Jie	G221HS218	60	from supplier
PhytoComm	Jing Jie	G221H1038822	60	from supplier
PhytoComm	Jing Jie	G221H1038421	40	from supplier
PhytoComm	Jing Jie	G221H1038822	60	from supplier
PhytoComm	Jing Jie	G221H1038221	40	from supplier
PhytoComm	Jing Jie	G221H1038621	60	from supplier
PhytoComm	Jing Jie	G221H1038621	60	from supplier
PhytoComm	Jing Jie	G221H1038621	60	from supplier
PhytoComm	Jing Jie	G221H1038121	40	from supplier
PhytoComm	Jing Jie	G221H1038621	60	from supplier
PhytoComm	Jing Jie	G221H1038421	40	from supplier
PhytoComm	Ku Shen	G233H0908321	75	from supplier
PhytoComm	Ku Shen	G233H0908321	40	from supplier
PhytoComm	Ku Shen	G233H0908822	60	from supplier
PhytoComm	Ku Shen	G233H0908822	60	from supplier
PhytoComm	Ku Shen	G233H0908521	60	from supplier
PhytoComm	Ku Shen	G233H0908521	60	from supplier
PhytoComm	Ku Shen	G233H0908321	40	from supplier
PhytoComm	Lai Fu Zi	G209H1220421	60	from supplier
PhytoComm	Lai Fu Zi	G209H1220321	40	from supplier
PhytoComm	Lai Fu Zi	G209H1220321	40	from supplier
PhytoComm	Lai Fu Zi	G209H1220821	60	from supplier
PhytoComm	Lai Fu Zi	G209H1220621	60	from supplier
PhytoComm	Lai Fu Zi	G209H1220321	75	from supplier
PhytoComm	Lai Fu Zi	G209H1220821	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Lai Fu Zi	G209H1220621	60	from supplier
PhytoComm	Lian Qiao	G105H1132621	60	from supplier
PhytoComm	Lian Qiao	G105H1132721	60	from supplier
PhytoComm	Lian Qiao	G105H1132421	40	from supplier
PhytoComm	Lian Qiao	G105H1132721	60	from supplier
PhytoComm	Lian Qiao	G105H1132621	60	from supplier
PhytoComm	Lian Qiao	G105H1132521	60	from supplier
PhytoComm	Lian Qiao	G105H1132821	60	from supplier
PhytoComm	Lian Qiao	G105H1132821	60	from supplier
PhytoComm	Lian Qiao	G105H1132221	40	from supplier
PhytoComm	Lian Qiao	G105HS252RW1	60	from supplier
PhytoComm	Lian Qiao	G105HS252RW1	60	from supplier
PhytoComm	Lian Qiao	G105H1132321	40	from supplier
PhytoComm	Lian Qiao	G105H1132421	75	from supplier
PhytoComm	Lian Qiao	G105H1132421	60	from supplier
PhytoComm	Lian Zi	G171H1501421	60	from supplier
PhytoComm	Lian Zi	G171HS312QM1	60	from supplier
PhytoComm	Lian Zi	G171HS312QM1	60	from supplier
PhytoComm	Lian Zi	G171HS312SK1	60	from supplier
PhytoComm	Lian Zi	G171HS312SK1	60	from supplier
PhytoComm	Lian Zi	G171H1501421	40	from supplier
PhytoComm	Lian Zi	G171H1501421	75	from supplier
PhytoComm	Ling Zhi	G109FB2405521	60	from supplier
PhytoComm	Ling Zhi	G109FB2405321	40	from supplier
PhytoComm	Ling Zhi	G109FB2405621	60	from supplier
PhytoComm	Ling Zhi	G109HS355RP1	60	from supplier
PhytoComm	Ling Zhi	G109FB2405621	60	from supplier
PhytoComm	Ling Zhi	G109FB2405321	40	from supplier
PhytoComm	Ling Zhi	G109HS355RP1	60	from supplier
PhytoComm	Ling Zhi	G109FB2405321	75	from supplier
PhytoComm	Long Dan (Cao)	G112H1705621	61	from supplier
PhytoComm	Long Dan (Cao)	G112H1705621	60	from supplier
PhytoComm	Long Dan (Cao)	G112H1705321	40	from supplier
PhytoComm	Long Dan (Cao)	G112H1705321	75	from supplier
PhytoComm	Long Dan (Cao)	G112H1705821	60	from supplier
PhytoComm	Long Dan (Cao)	G112H1705321	40	from supplier
PhytoComm	Long Dan (Cao)	G112H1705821	60	from supplier
PhytoComm	Long Yan Rou	G144H1708421	40	from supplier
PhytoComm	Long Yan Rou	G144H1708221	40	from supplier
PhytoComm	Long Yan Rou	G144HS324PW1	60	from supplier
PhytoComm	Long Yan Rou	G144H1708221	75	from supplier
PhytoComm	Long Yan Rou	G144HS324PW1	60	from supplier
PhytoComm	Long Yan Rou	G144H1708421	60	from supplier
PhytoComm	Long Yan Rou	G144H1708421	75	from supplier
PhytoComm	Long Yan Rou	G144H1708221	20	from supplier
PhytoComm	Long Yan Rou	G144HS324RN1	60	from supplier
PhytoComm	Long Yan Rou	G144HS324RN1	60	from supplier
PhytoComm	Long Yan Rou	G144H1708421	60	from supplier
PhytoComm	Lu Gen	G189H2014221	40	from supplier
PhytoComm	Lu Gen	G189H2014221	75	from supplier
PhytoComm	Lu Gen	G189H2014521	60	from supplier
PhytoComm	Lu Gen	G189H2014821	60	from supplier
PhytoComm	Lu Gen	G189H2014221	40	from supplier
PhytoComm	Lu Gen	G189H2014821	60	from supplier
PhytoComm	Ma Bian Cao	G248H1085823	60	from supplier
PhytoComm	Ma Bian Cao	G248H1085621	60	from supplier
PhytoComm	Ma Bian Cao	G248H1085421	40	from supplier
PhytoComm	Ma Bian Cao	G248H1085421	40	from supplier
PhytoComm	Ma Bian Cao	G248H1085823	60	from supplier
PhytoComm	Ma Bian Cao	G248H1085421	75	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Ma Huang	G097H1101623	60	from supplier
PhytoComm	Ma Huang	G097HS255SL1	60	from supplier
PhytoComm	Ma Huang	G097H1101422	60	from supplier
PhytoComm	Ma Huang	G097HS255SL1	60	from supplier
PhytoComm	Ma Huang	G097H1101522	60	from supplier
PhytoComm	Ma Huang	G097H1101422	75	from supplier
PhytoComm	Ma Huang	G097H1101623	60	from supplier
PhytoComm	Ma Huang	G097H1101223	30	from supplier
PhytoComm	Ma Huang	G097HS255RT1	60	from supplier
PhytoComm	Ma Huang	G097HS255RT1	60	from supplier
PhytoComm	Ma Huang	G097H1101422	40	from supplier
PhytoComm	Ma Huang Gen	G278H1102323	75	from supplier
PhytoComm	Ma Huang Gen	G278HS3930V1	60	from supplier
PhytoComm	Ma Huang Gen	G278H1102323	40	from supplier
PhytoComm	Ma Huang Gen	G278HS393RN1	60	from supplier
PhytoComm	Ma Huang Gen	G278H1102323	40	from supplier
PhytoComm	Ma Huang Gen	G278HS393RN1	60	from supplier
PhytoComm	Mai Men Dong	G176H1118221	40	from supplier
PhytoComm	Mai Men Dong	G176H1118521	60	from supplier
PhytoComm	Mai Men Dong	G176H1118621	60	from supplier
PhytoComm	Mai Men Dong	G176H1118621	60	from supplier
PhytoComm	Mai Men Dong	G176HS236SH1	60	from supplier
PhytoComm	Mai Men Dong	G176HS236SH1	60	from supplier
PhytoComm	Mai Men Dong	G176H1118622	60	from supplier
PhytoComm	Mai Men Dong	G176H1118521	60	from supplier
PhytoComm	Mai Men Dong	G176H1118322	40	from supplier
PhytoComm	Mai Men Dong	G176H1118622	60	from supplier
PhytoComm	Mai Ya	G123H1119821	60	from supplier
PhytoComm	Mai Ya	G123H1119421	40	from supplier
PhytoComm	Mai Ya	G123H1119521	60	from supplier
PhytoComm	Mai Ya	G123H1119521	60	from supplier
PhytoComm	Mai Ya	G123H1119521	60	from supplier
PhytoComm	Mai Ya	G123H1119421	75	from supplier
PhytoComm	Mai Ya	G123H1119321	40	from supplier
PhytoComm	Mai Ya	G123H1119821	60	from supplier
PhytoComm	Man Jing Zi	G250H1505622	60	from supplier
PhytoComm	Man Jing Zi	G250H1505622	60	from supplier
PhytoComm	Man Jing Zi	G250H1505822	60	from supplier
PhytoComm	Man Jing Zi	G250H1505522	60	from supplier
PhytoComm	Man Jing Zi	G250H1505121	30	from supplier
PhytoComm	Man Jing Zi	G250H1505822	60	from supplier
PhytoComm	Man Jing Zi	G250H1505421	60	from supplier
PhytoComm	Man Jing Zi	G250H1505421	40	from supplier
PhytoComm	Man Jing Zi	G250H1505421	75	from supplier
PhytoComm	Mang Xiao	G160H0627821	60	from supplier
PhytoComm	Mang Xiao	G160H0627421	75	from supplier
PhytoComm	Mang Xiao	G160H0627421	40	from supplier
PhytoComm	Mang Xiao	G160H0627821	60	from supplier
PhytoComm	Mang Xiao	G160H0627222	40	from supplier
PhytoComm	Mang Xiao	G160H0627623	60	from supplier
PhytoComm	Mang Xiao	G160H0627623	60	from supplier
PhytoComm	Mao Dong Qing	G319HS063RN1	60	from supplier
PhytoComm	Mao Dong Qing	G319HS063RN1	60	from supplier
PhytoComm	Mao Dong Qing	G319H0411221	40	from supplier
PhytoComm	Mao Dong Qing	G319H0411321	75	from supplier
PhytoComm	Mao Dong Qing	G319HS063SK1	60	from supplier
PhytoComm	Mao Dong Qing	G319HS063SK1	60	from supplier
PhytoComm	Mao Dong Qing	G319H0411321	40	from supplier
PhytoComm	Mao Dong Qing	G319H0411121	40	from supplier
PhytoComm	Mao Dong Qing	G319H0411321	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Mo Yao	G169H0732621	60	from supplier
PhytoComm	Mo Yao	G169H0732222	40	from supplier
PhytoComm	Mo Yao	G169H0732621	60	from supplier
PhytoComm	Mo Yao	G169H0732222	75	from supplier
PhytoComm	Mo Yao	G169H0732621	60	from supplier
PhytoComm	Mo Yao	G169H0732823	60	from supplier
PhytoComm	Mo Yao	G169H0732222	40	from supplier
PhytoComm	Mo Yao	G169H0732823	60	from supplier
PhytoComm	Mo Yao	G169H0732521	60	from supplier
PhytoComm	Mo Yao	G169H0732621	60	from supplier
PhytoComm	Mu Dan Pi	G166H0741621	60	from supplier
PhytoComm	Mu Dan Pi	G166H0741621	60	from supplier
PhytoComm	Mu Dan Pi	G166H0741421	60	from supplier
PhytoComm	Mu Dan Pi	G166HS129RT1	60	from supplier
PhytoComm	Mu Dan Pi	G166HS129RT1	60	from supplier
PhytoComm	Mu Dan Pi	G166H0741122	40	from supplier
PhytoComm	Mu Dan Pi	G166H0741321	40	from supplier
PhytoComm	Mu Dan Pi	G166H0741821	60	from supplier
PhytoComm	Mu Dan Pi	G166H0741321	75	from supplier
PhytoComm	Mu Dan Pi	G166H0741821	60	from supplier
PhytoComm	Mu Dan Pi	G166H0741321	40	from supplier
PhytoComm	Mu Gua	G063H0424521	60	from supplier
PhytoComm	Mu Gua	G063H0424321	60	from supplier
PhytoComm	Mu Gua	G063H0424221	75	from supplier
PhytoComm	Mu Gua	G063H0424221	40	from supplier
PhytoComm	Mu Gua	G063H0424621	60	from supplier
PhytoComm	Mu Gua	G063H0424621	60	from supplier
PhytoComm	Nü Zhen Zi	G138H0345821	60	from supplier
PhytoComm	Nü Zhen Zi	G138H0345421	40	from supplier
PhytoComm	Nü Zhen Zi	G138H0345721	60	from supplier
PhytoComm	Nü Zhen Zi	G138H0345421	75	from supplier
PhytoComm	Nü Zhen Zi	G138H0345721	60	from supplier
PhytoComm	Nü Zhen Zi	G138HS010SK1	60	from supplier
PhytoComm	Nü Zhen Zi	G138H0345821	60	from supplier
PhytoComm	Nü Zhen Zi	G138HS010SK1	60	from supplier
PhytoComm	Nü Zhen Zi	G138H0345321	40	from supplier
PhytoComm	Pi Pa Ye	G099H0832221	40	from supplier
PhytoComm	Pi Pa Ye	G099HS166RP1	60	from supplier
PhytoComm	Pi Pa Ye	G099H0832521	60	from supplier
PhytoComm	Pi Pa Ye	G099HS166RP1	60	from supplier
PhytoComm	Pi Pa Ye	G099H0832221	75	from supplier
PhytoComm	Pu Gong Ying	G239H1315623	60	from supplier
PhytoComm	Pu Gong Ying	G239H1315623	60	from supplier
PhytoComm	Pu Gong Ying	G239H1315623	60	from supplier
PhytoComm	Pu Gong Ying	G239H1315623	60	from supplier
PhytoComm	Pu Gong Ying	G239HS308SH1	60	from supplier
PhytoComm	Pu Gong Ying	G239HS308RN1	60	from supplier
PhytoComm	Pu Gong Ying	G239HS308RN1	60	from supplier
PhytoComm	Pu Gong Ying	G239H1315321	75	from supplier
PhytoComm	Pu Gong Ying	G239H1315321	40	from supplier
PhytoComm	Pu Gong Ying	G239H1315321	40	from supplier
PhytoComm	Pu Gong Ying	G239HS308SH1	60	from supplier
PhytoComm	Qiang Huo	G172H0840321	40	from supplier
PhytoComm	Qiang Huo	G172H0840521	60	from supplier
PhytoComm	Qiang Huo	G172H0840521	60	from supplier
PhytoComm	Qiang Huo	G172HS154SL1	60	from supplier
PhytoComm	Qiang Huo	G172H0840421	60	from supplier
PhytoComm	Qiang Huo	G172H0840422	60	from supplier
PhytoComm	Qiang Huo	G172HS154SL1	60	from supplier
PhytoComm	Qiang Huo	G172H0840621	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Qiang Huo	G172H0840721	60	from supplier
PhytoComm	Qiang Huo	G172H0840721	60	from supplier
PhytoComm	Qiang Huo	G172H0840621	60	from supplier
PhytoComm	Qin Jiao	G113H1059321	40	from supplier
PhytoComm	Qin Jiao	G113H1059521	60	from supplier
PhytoComm	Qin Jiao	G113H1059321	40	from supplier
PhytoComm	Qin Jiao	G113H1059321	75	from supplier
PhytoComm	Qin Jiao	G113H1059521	60	from supplier
PhytoComm	Qin Jiao	G113H1059321	60	from supplier
PhytoComm	Qing Hao	G031H0810321	40	from supplier
PhytoComm	Qing Hao	G031H0810721	60	from supplier
PhytoComm	Qing Hao	G031H0810721	60	from supplier
PhytoComm	Qing Hao	G031H0810321	75	from supplier
PhytoComm	Qing Hao	G031H0810321	40	from supplier
PhytoComm	Qing Pi	G074H0809323	40	from supplier
PhytoComm	Qing Pi	G074H0809323	40	from supplier
PhytoComm	Qing Pi	G074H0809621	60	from supplier
PhytoComm	Qing Pi	G074H0809821	60	from supplier
PhytoComm	Qing Pi	G074H0809821	60	from supplier
PhytoComm	Qing Pi	G074H0809323	75	from supplier
PhytoComm	Qing Pi	G074H0809621	60	from supplier
PhytoComm	Qing Pi	G074HS361SK1	60	from supplier
PhytoComm	Qing Pi	G074HS361SK1	60	from supplier
PhytoComm	Ren Dong Teng	G297H0751222	40	from supplier
PhytoComm	Ren Dong Teng	G297HS144SK1	60	from supplier
PhytoComm	Ren Dong Teng	G297HS144SK1	60	from supplier
PhytoComm	Rou Cong Rong	G072H1460621	60	from supplier
PhytoComm	Rou Cong Rong	G072HS108SG1	60	from supplier
PhytoComm	Rou Cong Rong	G072H1460121	20	from supplier
PhytoComm	Rou Cong Rong	G072H1460321	44	from supplier
PhytoComm	Rou Cong Rong	G072HS108SG1	60	from supplier
PhytoComm	Rou Cong Rong	G072H1460321	75	from supplier
PhytoComm	Rou Cong Rong	G072H1460321	40	from supplier
PhytoComm	Rou Gui	G069HS204RM1	60	from supplier
PhytoComm	Rou Gui	G069HS204RM1	60	from supplier
PhytoComm	Rou Gui	G069HS204QM1	60	from supplier
PhytoComm	Rou Gui	G069HS204QM1	60	from supplier
PhytoComm	Rou Gui	G069H0617123	75	from supplier
PhytoComm	Rou Gui	G069HS204RN1	60	from supplier
PhytoComm	Rou Gui	G069H0617121	20	from supplier
PhytoComm	Rou Gui	G069H0617123	40	from supplier
PhytoComm	Rou Gui	G069H0617323	60	from supplier
PhytoComm	Rou Gui	G069HS204RN1	60	from supplier
PhytoComm	Ru Xiang	G174H0820421	40	from supplier
PhytoComm	Ru Xiang	G174H0820421	40	from supplier
PhytoComm	Ru Xiang	G174H0820321	40	from supplier
PhytoComm	Ru Xiang	G174H0820421	60	from supplier
PhytoComm	Ru Xiang	G174H0820621	60	from supplier
PhytoComm	Ru Xiang	G174H0820421	75	from supplier
PhytoComm	Ru Xiang	G174H0820621	60	from supplier
PhytoComm	San Leng	G270H0311721	60	from supplier
PhytoComm	San Leng	G270H0311522	60	from supplier
PhytoComm	San Leng	G270H0311721	60	from supplier
PhytoComm	San Leng	G270H0311522	60	from supplier
PhytoComm	San Leng	G270H0311221	40	from supplier
PhytoComm	San Leng	G270H0311521	60	from supplier
PhytoComm	San Qi	G203H0312621	60	from supplier
PhytoComm	San Qi	G203H0312421	60	from supplier
PhytoComm	San Qi	G203HS030SL1	60	from supplier
PhytoComm	San Qi	G203H0312421	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	San Qi	G203H0312721	60	from supplier
PhytoComm	San Qi	G203H0312321	40	from supplier
PhytoComm	San Qi	G203HS030SL1	60	from supplier
PhytoComm	San Qi	G203H0312721	60	from supplier
PhytoComm	San Qi	G203H0312321	75	from supplier
PhytoComm	San Qi	G203H0312621	60	from supplier
PhytoComm	San Qi	G203H0312321	40	from supplier
PhytoComm	San Qi	G203H0312321	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042721	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042521	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042721	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042322	40	from supplier
PhytoComm	Sang Bai Pi	G161H1042322	40	from supplier
PhytoComm	Sang Bai Pi	G161H1042822	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042822	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042621	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042821	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042821	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042321	40	from supplier
PhytoComm	Sang Bai Pi	G161H1042621	60	from supplier
PhytoComm	Sang Bai Pi	G161H1042322	75	from supplier
PhytoComm	Sang Ji Shend	G146H1045822	60	from supplier
PhytoComm	Sang Ji Shend	G146H1045822	60	from supplier
PhytoComm	Sang Ji Shend	G146H1045621	60	from supplier
PhytoComm	Sang Ji Shend	G146H1045621	60	from supplier
PhytoComm	Sang Ji Shend	G146H1045821	60	from supplier
PhytoComm	Sang Ji Shend	G146H1045422	40	from supplier
PhytoComm	Sang Ji Shend	G146H1045422	75	from supplier
PhytoComm	Sang Ji Shend	G146H1045222	40	from supplier
PhytoComm	Sang Ji Shend	G146H1045322	40	from supplier
PhytoComm	Sang Ji Shend	G146H1045522	60	from supplier
PhytoComm	Sang Ji Shend	G146H1045821	60	from supplier
PhytoComm	Sang Ye	G162H1043421	75	from supplier
PhytoComm	Sang Ye	G162HS206RQ1	60	from supplier
PhytoComm	Sang Ye	G162H1043621	60	from supplier
PhytoComm	Sang Ye	G162H1043621	60	from supplier
PhytoComm	Sang Ye	G162HS206RQ1	60	from supplier
PhytoComm	Sang Ye	G162H1043421	40	from supplier
PhytoComm	Sang Zhi	G164H1044621	60	from supplier
PhytoComm	Sang Zhi	G164H1044222	40	from supplier
PhytoComm	Sang Zhi	G164H1044821	60	from supplier
PhytoComm	Sang Zhi	G164H1044222	40	from supplier
PhytoComm	Sang Zhi	G164H1044821	60	from supplier
PhytoComm	Sang Zhi	G164H1044721	60	from supplier
PhytoComm	Sang Zhi	G164H1044721	60	from supplier
PhytoComm	Sang Zhi	G164H1044222	75	from supplier
PhytoComm	Sang Zhi	G164H1044621	60	from supplier
PhytoComm	Sang Zhi	G164H1044421	60	from supplier
PhytoComm	Sha Ren	G018H0939321	75	from supplier
PhytoComm	Sha Ren	G018H0939821	60	from supplier
PhytoComm	Sha Ren	G018H0939821	60	from supplier
PhytoComm	Sha Ren	G018H0939622	60	from supplier
PhytoComm	Sha Ren	G018H0939522	60	from supplier
PhytoComm	Sha Ren	G018H0939622	60	from supplier
PhytoComm	Sha Ren	G018H0939322	40	from supplier
PhytoComm	Sha Ren	G018H0939321	40	from supplier
PhytoComm	Sha Ren	G018H0939321	40	from supplier
PhytoComm	Sha Ren	G018H0939523	60	from supplier
PhytoComm	Sha Ren	G018H0939523	60	from supplier
PhytoComm	Sha Shen (Bei)	G117H0754821	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Sha Shen (Bei)	G117H0754821	59	from supplier
PhytoComm	Sha Shen (Bei)	G117H0754322	40	from supplier
PhytoComm	Sha Shen (Bei)	G117H0754322	75	from supplier
PhytoComm	Shan Yao	G091H0330822	60	from supplier
PhytoComm	Shan Yao	G091H0330421	40	from supplier
PhytoComm	Shan Yao	G091H0330622	60	from supplier
PhytoComm	Shan Yao	G091H420302204	40	from supplier
PhytoComm	Shan Yao	G091H0330821	60	from supplier
PhytoComm	Shan Yao	G091H0330821	60	from supplier
PhytoComm	Shan Yao	G091HS024RQ1	60	from supplier
PhytoComm	Shan Yao	G091H420302204	75	from supplier
PhytoComm	Shan Yao	G091HS024RQ1	60	from supplier
PhytoComm	Shan Yao	G091H0330421	75	from supplier
PhytoComm	Shan Yao	G091H0330822	60	from supplier
PhytoComm	Shan Yao	G091H0330521	60	from supplier
PhytoComm	Shan Yao	G091H420302204	40	from supplier
PhytoComm	Shan Yao	G091H0330421	40	from supplier
PhytoComm	Shan Yao	G091H0330622	60	from supplier
PhytoComm	Shan Yu Rou	G080H0333421	60	from supplier
PhytoComm	Shan Yu Rou	G080HS027SK1	60	from supplier
PhytoComm	Shan Yu Rou	G080H0333321	40	from supplier
PhytoComm	Shan Yu Rou	G080H0333221	40	from supplier
PhytoComm	Shan Yu Rou	G080H0333621	60	from supplier
PhytoComm	Shan Yu Rou	G080H0333321	40	from supplier
PhytoComm	Shan Yu Rou	G080HS027RH1	60	from supplier
PhytoComm	Shan Yu Rou	G080HS027RH1	60	from supplier
PhytoComm	Shan Yu Rou	G080H0333321	75	from supplier
PhytoComm	Shan Yu Rou	G080H0333621	60	from supplier
PhytoComm	Shan Yu Rou	G080HS027SK1	60	from supplier
PhytoComm	She Gan	G048H1034422	60	from supplier
PhytoComm	She Gan	G048H1034621	60	from supplier
PhytoComm	She Gan	G048H1034621	60	from supplier
PhytoComm	She Gan	G048H1034221	75	from supplier
PhytoComm	She Gan	G048H1034221	40	from supplier
PhytoComm	She Gan	G048H1034221	20	from supplier
PhytoComm	Shen Qu	G156H0935421	40	from supplier
PhytoComm	Shen Qu	G156H0935421	40	from supplier
PhytoComm	Shen Qu	G156H0935621	60	from supplier
PhytoComm	Shen Qu	G156H0935821	60	from supplier
PhytoComm	Shen Qu	G156H0935821	60	from supplier
PhytoComm	Shen Qu	G156H0935421	75	from supplier
PhytoComm	Shen Qu	G156H0935621	60	from supplier
PhytoComm	Shen Qu	G156H0935521	60	from supplier
PhytoComm	Shen Qu	G156HS222RP1	60	from supplier
PhytoComm	Shen Qu	G156HS222RP1	60	from supplier
PhytoComm	Shen Qu	G156H0935522	60	from supplier
PhytoComm	Shen Qu	G156H0935522	60	from supplier
PhytoComm	Shen Qu	G156H0935121	40	from supplier
PhytoComm	Sheng Jiang	G253H0534221	40	from supplier
PhytoComm	Sheng Jiang	G253H0534623	60	from supplier
PhytoComm	Sheng Jiang	G253H0534421	40	from supplier
PhytoComm	Sheng Jiang	G253HS085SK1	60	from supplier
PhytoComm	Sheng Jiang	G253H0534623	60	from supplier
PhytoComm	Sheng Jiang	G253H0534522	60	from supplier
PhytoComm	Sheng Jiang	G253H0534421	40	from supplier
PhytoComm	Sheng Jiang	G253HS085SK1	60	from supplier
PhytoComm	Sheng Jiang	G253H0534321	40	from supplier
PhytoComm	Sheng Jiang	G253H0534421	75	from supplier
PhytoComm	Sheng Ma	G068H0437621	60	from supplier
PhytoComm	Sheng Ma	G068H0437621	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Sheng Ma	G068H0437421	40	from supplier
PhytoComm	Sheng Ma	G068H0437421	75	from supplier
PhytoComm	Shi Gao	G120H0541522	60	from supplier
PhytoComm	Shi Gao	G120H0541421	75	from supplier
PhytoComm	Shi Gao	G120H0541223	40	from supplier
PhytoComm	Shi Gao	G120H0541421	40	from supplier
PhytoComm	Shi Gao	G120H0541322	20	from supplier
PhytoComm	Shi Gao	G120HS081SH1	60	from supplier
PhytoComm	Shi Gao	G120H0541721	45	from supplier
PhytoComm	Shi Gao	G120H0541721	45	from supplier
PhytoComm	Shi Gao	G120HS081SH1	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601822	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601521	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601421	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601822	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601421	40	from supplier
PhytoComm	Shu Di (Huang)	G210H1601521	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601721	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601823	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601823	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601421	75	from supplier
PhytoComm	Shu Di (Huang)	G210H1601521	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601721	60	from supplier
PhytoComm	Shu Di (Huang)	G210H1601221	40	from supplier
PhytoComm	Suan Zao Ren	G254HS310SM1	60	from supplier
PhytoComm	Suan Zao Ren	G254HS310SM1	60	from supplier
PhytoComm	Suan Zao Ren	G254H1418121	10	from supplier
PhytoComm	Suan Zao Ren	G254H1418621	60	from supplier
PhytoComm	Suan Zao Ren	G254H1418621	60	from supplier
PhytoComm	Suan Zao Ren	G254H1418421	75	from supplier
PhytoComm	Suan Zao Ren	G254H1418622	60	from supplier
PhytoComm	Suan Zao Ren	G254H1418622	60	from supplier
PhytoComm	Suan Zao Ren	G254H1418421	60	from supplier
PhytoComm	Suan Zao Ren	G254H1418421	40	from supplier
PhytoComm	Suan Zao Ren	G254H1418322	40	from supplier
PhytoComm	Tai Zi Shen	G204HM066SK1	60	from supplier
PhytoComm	Tai Zi Shen	G204H0446421	40	from supplier
PhytoComm	Tai Zi Shen	G204H0446221	40	from supplier
PhytoComm	Tai Zi Shen	G204HM066SK1	60	from supplier
PhytoComm	Tai Zi Shen	G204H0446421	75	from supplier
PhytoComm	Tai Zi Shen	G204H0446721	60	from supplier
PhytoComm	Tai Zi Shen	G204H0446421	40	from supplier
PhytoComm	Tai Zi Shen	G204H0446721	60	from supplier
PhytoComm	Tai Zi Shen	G204H0446621	60	from supplier
PhytoComm	Tai Zi Shen	G204H0446621	60	from supplier
PhytoComm	Tai Zi Shen	G204H0446521	60	from supplier
PhytoComm	Tao Ren	G184HS220	60	from supplier
PhytoComm	Tao Ren	G184H1007321	40	from supplier
PhytoComm	Tao Ren	G184H1007321	75	from supplier
PhytoComm	Tao Ren	G184H1007321	30	from supplier
PhytoComm	Tao Ren	G184HS220SK1	60	from supplier
PhytoComm	Tao Ren	G184H1007223	40	from supplier
PhytoComm	Tao Ren	G184H1007621	60	from supplier
PhytoComm	Tao Ren	G184H1007621	60	from supplier
PhytoComm	Tao Ren	G184HS220SK1	60	from supplier
PhytoComm	Tao Ren	G184H1007422	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417621	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417821	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417821	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417621	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Tian Hua Fen	G242H0417421	75	from supplier
PhytoComm	Tian Hua Fen	G242H0417423	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417522	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417421	40	from supplier
PhytoComm	Tian Hua Fen	G242HS040SG1	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417421	40	from supplier
PhytoComm	Tian Hua Fen	G242HS040SG1	60	from supplier
PhytoComm	Tian Hua Fen	G242H0417321	50	from supplier
PhytoComm	Tian Men Dong	G037H0415321	40	from supplier
PhytoComm	Tian Men Dong	G037H0415221	40	from supplier
PhytoComm	Tian Men Dong	G037HS041RR1	60	from supplier
PhytoComm	Tian Men Dong	G037HS041RR1	60	from supplier
PhytoComm	Tian Men Dong	G037H0415321	40	from supplier
PhytoComm	Tian Men Dong	G037H0415121	40	from supplier
PhytoComm	Tian Men Dong	G037H0415521	60	from supplier
PhytoComm	Tian Men Dong	G037H0415321	75	from supplier
PhytoComm	Tu Fu Ling	G231H0325321	40	from supplier
PhytoComm	Tu Fu Ling	G231H0325621	60	from supplier
PhytoComm	Tu Fu Ling	G231HS009SK1	60	from supplier
PhytoComm	Tu Fu Ling	G231HS009RK1	60	from supplier
PhytoComm	Tu Fu Ling	G231H0325121	40	from supplier
PhytoComm	Tu Fu Ling	G231HS009SK1	60	from supplier
PhytoComm	Tu Fu Ling	G231H0325421	75	from supplier
PhytoComm	Tu Fu Ling	G231H0325421	60	from supplier
PhytoComm	Tu Fu Ling	G231HS009RK1	60	from supplier
PhytoComm	Tu Fu Ling	G231H0325621	60	from supplier
PhytoComm	Tu Fu Ling	G231H0325421	40	from supplier
PhytoComm	Tu Si Zi	G086H1243322	75	from supplier
PhytoComm	Tu Si Zi	G086H1243322	40	from supplier
PhytoComm	Tu Si Zi	G086H1243521	60	from supplier
PhytoComm	Tu Si Zi	G086H1243823	60	from supplier
PhytoComm	Tu Si Zi	G086H1243322	40	from supplier
PhytoComm	Tu Si Zi	G086H1243823	60	from supplier
PhytoComm	Tu Si Zi	G086H1243422	40	from supplier
PhytoComm	Tu Si Zi	G086H1243422	75	from supplier
PhytoComm	Tu Si Zi	G086H1243422	40	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456421	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281HS057QM1	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281HS057QM1	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456521	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456421	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456421	40	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456421	75	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456821	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456821	60	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456322	39	from supplier
PhytoComm	Wang Bu Liu Xing	G281H0456221	40	from supplier
PhytoComm	Wei Ling Xian	G075H0943422	75	from supplier
PhytoComm	Wei Ling Xian	G075H0943223	75	from supplier
PhytoComm	Wei Ling Xian	G075H0943521	60	from supplier
PhytoComm	Wei Ling Xian	G075H0943521	60	from supplier
PhytoComm	Wei Ling Xian	G075H0943223	40	from supplier
PhytoComm	Wei Ling Xian	G075H0943822	60	from supplier
PhytoComm	Wei Ling Xian	G075H0943822	60	from supplier
PhytoComm	Wei Ling Xian	G075H0943422	60	from supplier
PhytoComm	Wei Ling Xian	G075H0943422	40	from supplier
PhytoComm	Wei Ling Xian	G075H0943223	20	from supplier
PhytoComm	Wu Mei	G167HS192SK1	60	from supplier
PhytoComm	Wu Mei	G167H1026321	75	from supplier
PhytoComm	Wu Mei	G167H1026321	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Wu Mei	G167H1026521	60	from supplier
PhytoComm	Wu Mei	G167H1026221	40	from supplier
PhytoComm	Wu Mei	G167H1026521	60	from supplier
PhytoComm	Wu Mei	G167HS192SK1	60	from supplier
PhytoComm	Wu Mei	G167H1026421	60	from supplier
PhytoComm	Wu Mei	G167H1026321	40	from supplier
PhytoComm	Wu Mei	G167H1026821	60	from supplier
PhytoComm	Wu Mei	G167H1026821	60	from supplier
PhytoComm	Wu Wei Zi	G220H0401421	75	from supplier
PhytoComm	Wu Wei Zi	G220HS048RP1	60	from supplier
PhytoComm	Wu Wei Zi	G220H0401621	60	from supplier
PhytoComm	Wu Wei Zi	G220H0401621	60	from supplier
PhytoComm	Wu Wei Zi	G220H0401421	40	from supplier
PhytoComm	Wu Wei Zi	G220HS048SN1	60	from supplier
PhytoComm	Wu Wei Zi	G220H0401321	40	from supplier
PhytoComm	Wu Wei Zi	G220HS048SK1	60	from supplier
PhytoComm	Wu Wei Zi	G220HS048RP1	60	from supplier
PhytoComm	Wu Wei Zi	G220H0401122	40	from supplier
PhytoComm	Wu Wei Zi	G220HS048SK1	60	from supplier
PhytoComm	Wu Wei Zi	G220HS048SN1	60	from supplier
PhytoComm	Wu Yao	G140H1025221	20	from supplier
PhytoComm	Wu Yao	G140H1025421	40	from supplier
PhytoComm	Wu Yao	G140H1025721	60	from supplier
PhytoComm	Wu Yao	G140H1025421	40	from supplier
PhytoComm	Wu Yao	G140H1025421	75	from supplier
PhytoComm	Wu Yao	G140H1025521	60	from supplier
PhytoComm	Wu Yao	G140H1025721	60	from supplier
PhytoComm	Xi Xian Cao	G229F110311	40	from supplier
PhytoComm	Xi Xian Cao	G229F110311	75	from supplier
PhytoComm	Xi Xian Cao	G229F110311	40	from supplier
PhytoComm	Xi Xian Cao	G229H1431121	40	from supplier
PhytoComm	Xi Xian Cao	G229HS309SK1	60	from supplier
PhytoComm	Xi Xian Cao	G229HS309SK1	60	from supplier
PhytoComm	Xi Xian Cao	G229H1431321	60	from supplier
PhytoComm	Xi Xian Cao	G229H1431521	60	from supplier
PhytoComm	Xi Xian Cao	G229H1431521	60	from supplier
PhytoComm	Xia Ku Cao	G202H1064221	40	from supplier
PhytoComm	Xia Ku Cao	G202HS232RM1	60	from supplier
PhytoComm	Xia Ku Cao	G202HS232RM1	60	from supplier
PhytoComm	Xia Ku Cao	G202H1064421	75	from supplier
PhytoComm	Xia Ku Cao	G202HS232QM1	60	from supplier
PhytoComm	Xia Ku Cao	G202H1064421	60	from supplier
PhytoComm	Xia Ku Cao	G202HS232QM1	60	from supplier
PhytoComm	Xia Ku Cao	G202H1064421	20	from supplier
PhytoComm	Xia Ku Cao	G202H1064421	40	from supplier
PhytoComm	Xian He Cao	G008H0547321	75	from supplier
PhytoComm	Xian He Cao	G008H0547921	60	from supplier
PhytoComm	Xian He Cao	G008H0547321	40	from supplier
PhytoComm	Xian He Cao	G008H0547321	40	from supplier
PhytoComm	Xian He Cao	G008H0547321	40	from supplier
PhytoComm	Xian He Cao	G008H0547921	60	from supplier
PhytoComm	Xiang Fu	G088H0930422	40	from supplier
PhytoComm	Xiang Fu	G088H0930821	60	from supplier
PhytoComm	Xiang Fu	G088H0930222	40	from supplier
PhytoComm	Xiang Fu	G088H0930621	60	from supplier
PhytoComm	Xiang Fu	G088H0930621	60	from supplier
PhytoComm	Xiang Fu	G088H0930422	75	from supplier
PhytoComm	Xiang Fu	G088H0930322	40	from supplier
PhytoComm	Xiang Fu	G088H0930521	60	from supplier
PhytoComm	Xiang Fu	G088H0930621	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Xiang Fu	G088H0930621	60	from supplier
PhytoComm	Xiang Fu	G088H0930321	40	from supplier
PhytoComm	Xiang Fu	G088H0930821	60	from supplier
PhytoComm	Xie Bai	G504H1830422	60	from supplier
PhytoComm	Xie Bai	G504H1830421	40	from supplier
PhytoComm	Xie Bai	G504H1830321	50	from supplier
PhytoComm	Xie Bai	G504H1830821	60	from supplier
PhytoComm	Xie Bai	G504H1830621	60	from supplier
PhytoComm	Xie Bai	G504H1830421	40	from supplier
PhytoComm	Xie Bai	G504H1830821	60	from supplier
PhytoComm	Xie Bai	G504H1830621	60	from supplier
PhytoComm	Xie Bai	G504H1830421	75	from supplier
PhytoComm	Xin Yi	G153H0758523	60	from supplier
PhytoComm	Xin Yi	G153H0758523	60	from supplier
PhytoComm	Xin Yi	G153H0758321	40	from supplier
PhytoComm	Xin Yi	G153HS139SH1	60	from supplier
PhytoComm	Xin Yi	G153H0758821	60	from supplier
PhytoComm	Xin Yi	G153H0758321	75	from supplier
PhytoComm	Xin Yi	G153H0758521	60	from supplier
PhytoComm	Xin Yi	G153H0758322	45	from supplier
PhytoComm	Xin Yi	G153H0758721	60	from supplier
PhytoComm	Xin Yi	G153HS139SH1	60	from supplier
PhytoComm	Xin Yi	G153H0758321	40	from supplier
PhytoComm	Xin Yi	G153H0758721	60	from supplier
PhytoComm	Xin Yi	G153H0758821	60	from supplier
PhytoComm	Xu Duan	G092H2106321	40	from supplier
PhytoComm	Xu Duan	G092H2106822	60	from supplier
PhytoComm	Xu Duan	G092H2106822	60	from supplier
PhytoComm	Xu Duan	G092H2106621	60	from supplier
PhytoComm	Xu Duan	G092H2106321	40	from supplier
PhytoComm	Xu Duan	G092H2106621	60	from supplier
PhytoComm	Xu Duan	G092H2106821	60	from supplier
PhytoComm	Xu Duan	G092H2106821	60	from supplier
PhytoComm	Xu Duan	G092H2106321	75	from supplier
PhytoComm	Xu Duan	G092H2106422	60	from supplier
PhytoComm	Xuan Shen	G223HS092QH1	60	from supplier
PhytoComm	Xuan Shen	G223HS092RN1	60	from supplier
PhytoComm	Xuan Shen	G223HS092RN1	60	from supplier
PhytoComm	Xuan Shen	G223H0537421	75	from supplier
PhytoComm	Xuan Shen	G223HS092SH1	60	from supplier
PhytoComm	Xuan Shen	G223HS092SH1	60	from supplier
PhytoComm	Xuan Shen	G223H0537521	60	from supplier
PhytoComm	Xuan Shen	G223H0537421	40	from supplier
PhytoComm	Xuan Shen	G223H0537521	60	from supplier
PhytoComm	Xuan Shen	G223HS092QH1	60	from supplier
PhytoComm	Xuan Shen	G223H0537322	40	from supplier
PhytoComm	Yan Hu Suo	G081H0805223	80	from supplier
PhytoComm	Yan Hu Suo	G081H0805422	40	from supplier
PhytoComm	Yan Hu Suo	G081H0805422	60	from supplier
PhytoComm	Yan Hu Suo	G081H0805821	60	from supplier
PhytoComm	Yan Hu Suo	G081HS161SG1	60	from supplier
PhytoComm	Yan Hu Suo	G081HS161SG1	60	from supplier
PhytoComm	Yan Hu Suo	G081H0805821	60	from supplier
PhytoComm	Yan Hu Suo	G081H0805521	60	from supplier
PhytoComm	Yan Hu Suo	G081H0805622	60	from supplier
PhytoComm	Yan Hu Suo	G081H0805622	60	from supplier
PhytoComm	Yan Hu Suo	G081H0805422	75	from supplier
PhytoComm	Ye Jiao Teng	G288H0885223	75	from supplier
PhytoComm	Ye Jiao Teng	G288H0885522	60	from supplier
PhytoComm	Ye Jiao Teng	G288H0885223	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Ye Jiao Teng	G288H0885621	60	from supplier
PhytoComm	Ye Jiao Teng	G288H0885522	60	from supplier
PhytoComm	Ye Jiao Teng	G288H0885621	60	from supplier
PhytoComm	Ye Jiao Teng	G288HS163RM1	60	from supplier
PhytoComm	Ye Jiao Teng	G288HS163RM1	60	from supplier
PhytoComm	Ye Jiao Teng	G288HS163SK1	60	from supplier
PhytoComm	Ye Jiao Teng	G288H0885223	40	from supplier
PhytoComm	Ye Jiao Teng	G288HS163SK1	60	from supplier
PhytoComm	Ye Jiao Teng	G288H0885522	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050421	40	from supplier
PhytoComm	Yi Mu Cao	G134H1050822	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050822	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050621	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050621	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050321	40	from supplier
PhytoComm	Yi Mu Cao	G134H1050521	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050421	75	from supplier
PhytoComm	Yi Mu Cao	G134H1050821	60	from supplier
PhytoComm	Yi Mu Cao	G134H1050821	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701522	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701222	40	from supplier
PhytoComm	Yi Yi Ren	G078HS336RL1	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701622	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701323	40	from supplier
PhytoComm	Yi Yi Ren	G078HS336RL1	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701622	60	from supplier
PhytoComm	Yi Yi Ren	G078HS336SH1	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701323	75	from supplier
PhytoComm	Yi Yi Ren	G078HS336SH1	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701522	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701421	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701522	60	from supplier
PhytoComm	Yi Yi Ren	G078H1701323	40	from supplier
PhytoComm	Yi Yi Ren	G078H1701522	60	from supplier
PhytoComm	Yin Chen Hao	G034H1030421	60	from supplier
PhytoComm	Yin Chen Hao	G034H1030121	40	from supplier
PhytoComm	Yin Chen Hao	G034H1030221	75	from supplier
PhytoComm	Yin Chen Hao	G034H1030221	40	from supplier
PhytoComm	Yin Chen Hao	G034H1030721	60	from supplier
PhytoComm	Yin Chen Hao	G034H1030621	60	from supplier
PhytoComm	Yin Chen Hao	G034H1030721	60	from supplier
PhytoComm	Yin Chen Hao	G034H1030221	40	from supplier
PhytoComm	Yin Yang Huo	G098H1114622	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114422	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114621	60	from supplier
PhytoComm	Yin Yang Huo	G098HS241SH1	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114621	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114622	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114821	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114223	75	from supplier
PhytoComm	Yin Yang Huo	G098H1114821	60	from supplier
PhytoComm	Yin Yang Huo	G098H1114223	40	from supplier
PhytoComm	Yin Yang Huo	G098HS241SH1	60	from supplier
PhytoComm	Yin Yang Huo	FA1114702	40	from supplier
PhytoComm	Yin Yang Huo	G098H1114223	40	from supplier
PhytoComm	Yu Jin	G084H2201821	60	from supplier
PhytoComm	Yu Jin	G084H2201221	20	from supplier
PhytoComm	Yu Jin	G084H2201821	60	from supplier
PhytoComm	Yu Jin	G084H2201521	60	from supplier
PhytoComm	Yu Jin	G084H2201621	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Yu Jin	G084H2201421	60	from supplier
PhytoComm	Yu Jin	G084H2201322	75	from supplier
PhytoComm	Yu Jin	G084H2201521	60	from supplier
PhytoComm	Yu Jin	G084H2201421	60	from supplier
PhytoComm	Yu Jin	G084H2201621	60	from supplier
PhytoComm	Yu Jin	G084HS357RN1	60	from supplier
PhytoComm	Yu Jin	G084HS357RN1	60	from supplier
PhytoComm	Yu Jin	G084H2201322	40	from supplier
PhytoComm	Yu Jin	G084H2201322	40	from supplier
PhytoComm	Yu Xing Cao	G124H1149823	60	from supplier
PhytoComm	Yu Xing Cao	G124H1149823	60	from supplier
PhytoComm	Yu Xing Cao	G124H1149222	40	from supplier
PhytoComm	Yu Xing Cao	G124H1149423	60	from supplier
PhytoComm	Yu Xing Cao	G124H1149222	40	from supplier
PhytoComm	Yu Xing Cao	G124H1149222	75	from supplier
PhytoComm	Yu Zhu	G195HS087SK1	60	from supplier
PhytoComm	Yu Zhu	G195H0517521	60	from supplier
PhytoComm	Yu Zhu	G195H0517521	60	from supplier
PhytoComm	Yu Zhu	G195H0517321	75	from supplier
PhytoComm	Yu Zhu	G195H0517321	40	from supplier
PhytoComm	Yu Zhu	G195H0517321	40	from supplier
PhytoComm	Yu Zhu	G195H0517421	60	from supplier
PhytoComm	Yu Zhu	G195H0517121	40	from supplier
PhytoComm	Yu Zhu	G195HS087SK1	60	from supplier
PhytoComm	Yuan Zhi	G194H1406522	60	from supplier
PhytoComm	Yuan Zhi	G194HS306SH1	60	from supplier
PhytoComm	Yuan Zhi	G194H1406622	60	from supplier
PhytoComm	Yuan Zhi	G194HS306SH1	60	from supplier
PhytoComm	Yuan Zhi	G194H1406622	60	from supplier
PhytoComm	Yuan Zhi	G194HS306RP1	60	from supplier
PhytoComm	Yuan Zhi	G194HS306RP1	60	from supplier
PhytoComm	Yuan Zhi	G194H1406222	40	from supplier
PhytoComm	Yuan Zhi	G194H1406222	40	from supplier
PhytoComm	Yuan Zhi	G194H1406222	75	from supplier
PhytoComm	Ze Xie	G011H1605122	40	from supplier
PhytoComm	Ze Xie	G011H1605422	60	from supplier
PhytoComm	Ze Xie	G011H1605822	60	from supplier
PhytoComm	Ze Xie	G011H1605822	60	from supplier
PhytoComm	Ze Xie	G011H1605422	75	from supplier
PhytoComm	Ze Xie	G011H1605821	60	from supplier
PhytoComm	Ze Xie	G011H1605422	60	from supplier
PhytoComm	Ze Xie	G011H1605422	40	from supplier
PhytoComm	Ze Xie	G011H1605821	60	from supplier
PhytoComm	Ze Xie	G011H1605521	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712321	40	from supplier
PhytoComm	Zhe Bei Mu	G108H0712321	75	from supplier
PhytoComm	Zhe Bei Mu	G108H0712721	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712422	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712721	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712321	40	from supplier
PhytoComm	Zhe Bei Mu	G108H0712821	60	from supplier
PhytoComm	Zhe Bei Mu	G108HS138SG1	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712821	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712121	40	from supplier
PhytoComm	Zhe Bei Mu	G108HS138SG1	60	from supplier
PhytoComm	Zhe Bei Mu	G108H0712521	60	from supplier
PhytoComm	Zhi Ke	G043H0921621	60	from supplier
PhytoComm	Zhi Ke	G043H0921221	75	from supplier
PhytoComm	Zhi Ke	G043H0921221	40	from supplier
PhytoComm	Zhi Ke	G043HS169SL1	60	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Zhi Ke	G043H0921422	60	from supplier
PhytoComm	Zhi Ke	G043HS169SL1	60	from supplier
PhytoComm	Zhi Ke	G043HS169RR1	60	from supplier
PhytoComm	Zhi Ke	G043H0921621	60	from supplier
PhytoComm	Zhi Ke	G043H0921221	40	from supplier
PhytoComm	Zhi Ke	G043HS169RR1	60	from supplier
PhytoComm	Zhi Mu	G019H0824521	60	from supplier
PhytoComm	Zhi Mu	G019HS153SG1	60	from supplier
PhytoComm	Zhi Mu	G019HS153SG1	60	from supplier
PhytoComm	Zhi Mu	G019H0824521	60	from supplier
PhytoComm	Zhi Mu	G019H0824821	60	from supplier
PhytoComm	Zhi Mu	G019H0824321	40	from supplier
PhytoComm	Zhi Mu	G019H0824821	60	from supplier
PhytoComm	Zhi Mu	G019H0824421	75	from supplier
PhytoComm	Zhi Mu	G019H0824322	40	from supplier
PhytoComm	Zhi Mu	G019HS153QL1	60	from supplier
PhytoComm	Zhi Mu	G019H0824421	40	from supplier
PhytoComm	Zhi Mu	G019HS153ON1	60	from supplier
PhytoComm	Zhi Mu	G019HS153QL1	60	from supplier
PhytoComm	Zhi Shi	G044H0922221	40	from supplier
PhytoComm	Zhi Shi	G044H0922821	60	from supplier
PhytoComm	Zhi Shi	G044HS168SK1	60	from supplier
PhytoComm	Zhi Shi	G044H0922821	60	from supplier
PhytoComm	Zhi Shi	G044H0922321	40	from supplier
PhytoComm	Zhi Shi	G044H0922321	75	from supplier
PhytoComm	Zhi Shi	G044H0922521	60	from supplier
PhytoComm	Zhi Shi	G044H0922522	60	from supplier
PhytoComm	Zhi Shi	G044H0922522	60	from supplier
PhytoComm	Zhi Shi	G044HS168SK1	60	from supplier
PhytoComm	Zhi Shi	G044H0922321	40	from supplier
PhytoComm	Zhi Zi	G110H1137321	60	from supplier
PhytoComm	Zhi Zi	G110H1137321	40	from supplier
PhytoComm	Zhi Zi	G110H1137321	40	from supplier
PhytoComm	Zhi Zi	G110H1137521	60	from supplier
PhytoComm	Zhi Zi	G110H1137321	60	from supplier
PhytoComm	Zhi Zi	G110H1137621	60	from supplier
PhytoComm	Zhi Zi	G110H1137523	60	from supplier
PhytoComm	Zhi Zi	G110H1137621	60	from supplier
PhytoComm	Zhi Zi	G110H1137523	60	from supplier
PhytoComm	Zhi Zi	G110H1137822	60	from supplier
PhytoComm	Zhi Zi	G110H1137822	60	from supplier
PhytoComm	Zhu Ru	G047H0625422	60	from supplier
PhytoComm	Zhu Ru	G047H0625122	20	from supplier
PhytoComm	Zhu Ru	G047H0625421	60	from supplier
PhytoComm	Zhu Ru	G047HS109QM2	60	from supplier
PhytoComm	Zhu Ru	G047HS109QM2	60	from supplier
PhytoComm	Zhu Ru	G047H0625421	75	from supplier
PhytoComm	Zhu Ru	G047H0625621	60	from supplier
PhytoComm	Zhu Ru	G047H0625621	60	from supplier
PhytoComm	Zhu Ru	G047H0625421	40	from supplier
PhytoComm	Zhu Ru	G047HS109SK1	60	from supplier
PhytoComm	Zhu Ru	G047H0625321	40	from supplier
PhytoComm	Zhu Ru	G047HS109SK1	60	from supplier
PhytoComm	Zhu Ru	G047H0625122	40	from supplier
PhytoComm	Zi Hua Di Ding	G249H1216621	60	from supplier
PhytoComm	Zi Hua Di Ding	G249H1216321	40	from supplier
PhytoComm	Zi Hua Di Ding	G249H1216621	60	from supplier
PhytoComm	Zi Hua Di Ding	G249H1216321	60	from supplier
PhytoComm	Zi Hua Di Ding	G249H1216321	75	from supplier
PhytoComm	Zi Hua Di Ding	G249H1216321	40	from supplier

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra	Certificate
PhytoComm	Zi Hua Di Ding	G249HS269SK1	60	from supplier
PhytoComm	Zi Hua Di Ding	G249HS269SK1	60	from supplier
PhytoComm	Zi Su Zi	G183H1213421	75	from supplier
PhytoComm	Zi Su Zi	G183H1213822	60	from supplier
PhytoComm	Zi Su Zi	G183H1213621	60	from supplier
PhytoComm	Zi Su Zi	G183H1213421	60	from supplier
PhytoComm	Zi Su Zi	G183H1213621	60	from supplier
PhytoComm	Zi Su Zi	G183H1213421	40	from supplier
PhytoComm	Zi Su Zi	G183H1213321	40	from supplier
PhytoComm	Zi Su Zi	G183H1213822	60	from supplier
PhytoComm	Zi Wan	G038H1215621	60	from supplier
PhytoComm	Zi Wan	G038H1215621	60	from supplier
PhytoComm	Zi Wan	G038H1215221	40	from supplier
PhytoComm	Zi Wan	G038H1215321	40	from supplier
PhytoComm	Zi Wan	G038H1215221	40	from supplier
PhytoComm	Zi Wan	G038H1215321	40	from supplier
PhytoComm	Zi Wan	G038H1215321	60	from supplier
PhytoComm	Zi Wan	G038H1215321	75	from supplier
PhytoComm	Zi Wan	G038H1215221	75	from supplier
Roth	Brilliant blue R...	314217301	60	from supplier
Roth	Deoxycholic acid...	385212633	60	not necessary
Roth	Deoxycholic acid...	154212633	60	not necessary
Roth	Deoxycholic acid...	084181622	40	1405618
Roth	Esculin	155206581	60	not necessary
Roth	Theobromine, 99%	145191010	60	AR-16-FG-002141-01
Roth	Theobromine, 99%	365191010	60	not necessary
Roth	Trypan blue (C.I...	204209355	60	from supplier
Safic-Alcan	Dimercaptosuccin...	20180130	40	from supplier
Sand 12 GmbH	D-ribose	L201211720	35	not necessary
Sigma-Aldrich	Bromelain	BCBM1589V	40	not necessary
Sigma-Aldrich	Bromothymol blue	MKBQ5404V	40	pending
Sigma-Aldrich	Dexamethasone is...	Y0000538	40	not necessary
Sigma-Aldrich	Dexamethasone so...	LRAC1690	40	not necessary
Sigma-Aldrich	Isoconazole	1.1	40	not necessary
Sigma-Aldrich	Isoconazole nitr...	1.1	40	not necessary
Sigma-Aldrich	Metronidazole be...	1.3	40	not necessary
Sigma-Aldrich	Prednisolone piv...	P2805000	40	not necessary
Sigma-Aldrich	Sertaconazole ni...	2.0	26	not necessary
Spinnrad	Inulin	1112300433	40	not necessary
Th. Geyer	Calcium colamine...	36986	40	from supplier
Th. Geyer	Mg-Ca-K colamine...	54610	40	from supplier
Th. Geyer	Mg-Ca-K colamine...	54610	40	from supplier
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	15-7221-80338	60	1316Q-07806
THC-Pharma GmbH	Cannabidiol	CBD-API17-003	40	not necessary
THC-Pharma GmbH	Cannabidiol	CBD-API17-003	31	not necessary
THC-Pharma GmbH	Cannabidiol	CBD-API16-003	40	2017B-03911
THC-Pharma GmbH	Cannabidiol	CBD-API16-004	40	2017A-03911
THC-Pharma GmbH	Cannabidiol	CBD-API1801	40	not necessary
THC-Pharma GmbH	Cannabidiol	CBD-API1801	40	not necessary
VMP Chemiekon...	D-(+)-glucosamin...	144212943	40	20140402
Yino Pharma	Heparin sodium	YN010121805010	40	not necessary
Yino Pharma	Heparin sodium	YN010111805002	40	not necessary
Yino Pharma	Heparin sodium	YN01011R1170803	40	not necessary
Yino Pharma	Heparin sodium	YN01011R1170803	40	not necessary
Yino Pharma	Sodium chondroit...	YN010511805005	40	not necessary
Yino Pharma	Sodium chondroit...	YN010511806007	40	not necessary

*The spectra were taken from a part of the same sample which the supplier used for its analysis (identity, purity, content). In this case the manufacturer certificates proof the identity and quality of the sample from those the spectra descend.

Appendix C: Additional validation samples (*Type C*)

Spectra for substances from the field which cannot be identified with this model also enter the validation. In this manner, it is verified that the model also rejects unknown substances. The spectra for these samples were recorded by *Apo-Ident* customers. They belong to *Type C*. The information provided by the customer regarding the manufacturer and batch number is taken over by *HiperScan GmbH* to a large extent unchecked.

The samples originate from 31 008 batches. From these, 94 638 spectra were recorded. The spectra recorded on independent samples of substances from the field which can be identified with the model are listed respectively in the section *Type C* for the individual substances and do not appear again elsewhere in this list.

Supplier	Substance	Batch	Spectra
/Sanacorp	Betamethasone valerate	17078701	14
/Sanacorp	Magnesium hydrogen ci...	90170	1
A	Oxytetracycline hydro...	130602202	2
A.H. 03.06.13 ek: 6,50EUR	Progesterone, microni...	1304014-01	1
A.H. 27.05.13 ek: 6,75EUR	Male speedwell wort, ...	274239	1
A.H.D	Minoxidil	1210056-01	21
A.H.D.	Emulsifying cetostear...	12L19-N03	9
A.H.D.	Lidocaine hydrochlori...	1202035-02	44
A.H.D.	Tannin	13E27-N08	2
ACM	Acetylsalicylic acid	KNC2077J14	1
ACM	Acetylsalicylic acid	KNC2152J13	1
ACM	Ascorbic acid	Knc1035j16	1
ACM	Ascorbic acid	KNC1035J16	2
ACM	Ascorbic acid	KNC1488J14	6
ACM	Benzocaine	KNC1643J15	1
ACM	Benzoic acid	KNC0558J15	1
ACM	Benzoic acid	KNC0396J17	1
ACM	Caffeine	KNC1294J14	2
ACM	Calcium citrate tetra...	KNC0740J15	1
ACM	Cocoa butter	KNC1605J15	1
ACM	Glucose monohydrate	KNC1164J15	1
ACM	Glucose monohydrate	KNC1282J14	1
ACM	Hydroquinone	KNC1047J15	2
ACM	Hydroquinone	KNC0635J17	1
ACM	Lidocaine hydrochlori...	KNC0625J17	1
ACM	Lidocaine hydrochlori...	KNC1247J15	1
ACM	Macrogols 4000	KNC0366J15	2
ACM	Magnesium hydrogen ci...	KNC0829J15	1
ACM	Magnesium hydrogen ci...	KNC1377J15	2
ACM	Metronidazole	EG:207-136-1	1
ACM	Potassium citrate	KNC0196J15	4
ACM	Salicylic acid	KNC0461J15	1
ACM	Urea	KNC0557J15	3
ACM	Urea	KNC0236J16	2
ACM	Urea	KNC1795J14	1
ACM / Herba	Macrogols 1500	KN-0826/10	1
ACM / Herba	Thymol	KNC1166J14	1
ACM / Phoenix	Theophylline-ethylene...	KNC0896J14	1
agron/Krieger	Betamethasone valerate	12B2B-N01	1
AHD	Bacitracin	17L11-B02-344244	1
AHD	Bitter fennel	15006502	1
AHD	Camphor	17288306	4
AHD	Cetomacrogol emulsify...	8073743	1
AHD	Chloramphenicol	18017807	4
AHD	Cholesterol	17K16-B11-352718	3
AHD	Clindamycin hydrochlo...	130402X5	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
AHD	Dexamethasone	17L19-B03-344099	1
AHD	Diclofenac sodium	15D14B04314450	1
AHD	Diltiazem hydrochloride	17043210	2
AHD	Diltiazem hydrochloride	18301108	3
AHD	Diltiazem hydrochloride	17337017	3
AHD	Erythromycin	18110119-01	1
AHD	Glucose monohydrate	14102407	3
AHD	Hard fat in pastilles	277212	3
AHD	Hydrocortisone	13069507	6
AHD	Lactose monohydrate	13I23-B03-291114	1
AHD	Lidocaine hydrochlori...	18179904	3
AHD	Lidocaine hydrochlori...	18179910	2
AHD	Macrogols 4000	17E02-B03	1
AHD	Macrogols 4000	15G15B07-313592	1
AHD	Macrogols 4000	17K06-B05-345076	3
AHD	Macrogols 4000	14L12-B06-301991	2
AHD	Macrogols 4000	17K06-B05	2
AHD	Macrogols 4000	11H11-B05340782	1
AHD	Macrogols 4000	18G30-B01-352394	2
AHD	Macrogols 4000	17E02-B03-335507	4
AHD	Menthol	18128014	1
AHD	Methoxsalen	16103019	1
AHD	Minoxidil	18149801	1
AHD	Minoxidil	15230812	1
AHD	Neomycin sulfate	14A31-B02-294698	4
AHD	Nystatin	1612001-01	1
AHD	Potassium sorbate	13065216	13
AHD	Prednisolone, microni...	160711-01	1
AHD	Prednisolone, microni...	121202103	2
AHD	Progesterone, microni...	13I11-N11	3
AHD	Propranolol hydrochlo...	941895	1
AHD	Propyl-4-hydroxybenzo...	11373523	3
AHD	Sodium benzoate	1711019-01	3
AHD	Sodium citrate	1710006-01	10
AHD	Sodium citrate	17323602	3
AHD	Sodium monohydrogenph...	112345013	1
AHD	Sucrose	18154602	2
AHD	Testosterone propionate	AP12E04-N01TP	1
AHD	Tetracaine hydrochlor...	18060614	1
AHD	Tragacanth, powdered	15379906	3
AHD	Tretinoin	16070002	3
AHD	Urea	3514E-02544	1
AHD	Wax, yellow	18299202	2
AHD 06.03.15 9,47	Triamcinolone acetonide	14355102	14
AHD 10.02.17 EK.3,50	Aniseed	16200801	4
AHD 15.06.16 EK:18,44	Betamethasone valerate	15398505	12
AHD 17.11.14 EK 7,60	Epinephrine hydrogen ...	Sana30.10.14EK260	1
AHD 20.06.14 EK 37,29	Mometasone furoate, m...	12C07-B05-291268	1
AHD 22.07.14	Salicylic acid	134002	1
AHD 22.12.16 EK 29,74	Triamcinolone acetonide	16D14B11322136	1
AHD 27.10 Ek.7,07	Clotrimazole	1607046	2
AHD 29.10.14 EK6,12	Gentamicin sulfate	13205004	1
ahd,11,95EUR 10.04.14	Hydrocortisone	13341504-07R07	1
AHD/Fargon	Betamethasone valerate	14H28-Bo1-302865	1
ahd/lamotte	Bismuth gallath, basic	13718_b13_292058	1
ahd/lamotte	Cetostearyl alcohol	13210504	1
AHD;11.05.15;3,23	Bitter fennel	14376408	3
AHDE 20.08.15	Chloramphenicol	DY00912128	3
AHL	Hydrocortisone	15125-B01-313084	1
AHL	Hydroquinone	16023102	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
AHL	Triclosan	1801049-01	10
AHZ	Clobetasol propionate	1606022	4
Akazie	Black cumin seed	20.09.2015	1
Akazie	Cinnamon, cut	20.09.2015	1
Akazie	Dill fruit	20.09.2015	1
Akazie	Linseed	20.09.2015	2
Akazie	Rosemary leaf, cut	20.09.2015	1
Akazien-Apo	Cinnamon, cut	Test	1
Aldrich	Sodium proprionate	U02796-056	1
Alfred Galke	Alchemilla, cut	20921	1
Alfred Galke	Bitter fennel	18270	1
Alfred Galke	Caraway fruit	20069	1
Alfred Galke	Melissa leaf, cut	D-NI-1-9516-BC	1
Alfred Galke	Raspberry leaf, cut	19851	1
Alfred Galke	Shepherd's purse, cut	20591	1
Alfred Galke	Yarrow herb, cut	19861	1
ALH	Chlortetracycline hyd...	16B10-B09-320304	2
ALH	Gentamicin sulfate	15258813	11
ALH	Metronidazole	1609028-01	35
ALH	Phenylephrine hydroch...	16A20-B07-322159	1
ALH	Sorbic acid	15382007	5
ALH	Tropicamide	15F18-B04-310971	1
alliance healthcare 5...	Estradiol benzoate	10118576	1
Allianc. Healthc. /c	Glucose	16050510	3
Alliance	Agrimony, cut	11374004	5
Alliance	Hawthorn flower	13010403	1
Alliance	Heather, cut	255389	1
Alliance	Melatonin	15254316	2
Alliance	Povidone 25, soluble	16307808	2
Alliance	Povidone 25, soluble	17209106	1
Alliance	Progesterone, microni...	15i11-b05-315224	1
Alliance	Progesterone, microni...	15I11-B05-315224	3
Alliance	Sildenafil citrate	31.07.18	1
Alliance	Sodium hydrogen carbo...	12324704	1
Alliance	Sodium monohydrogenph...	12256606	1
Alliance H.	Miconazole nitrate	5137910	1
Alliance health	Fructose	14K17-B99-300913	3
Alliance Health	Aluminium potassium s...	13025401	3
Alliance Health	Erythromycin	120503502	2
Alliance Health	Tetracaine hydrochlor...	14C18-B03-293967	2
Alliance Health 03.06...	Citric acid	13K18-B34	8
Alliance Health 05.05...	Tetracaine hydrochlor...	13H21-B02-286559	2
Alliance Health 11.01...	Hydrocortisone	1315505	1
Alliance Health 13.02...	Fructose	13C25-N03	1
Alliance Health 17.01...	Clobetasol propionate	APB00547/112010C	1
Alliance Health 22.08...	Tetracycline hydrochl...	13201307	4
Alliance Health 27.12...	Hypromellose 4000	13D03-N03	1
alliance health. 05.0...	Emulsifying cetostear...	14G31-B03-297610	4
Alliance Health25.04...	Hydrocortisone	13341503	4
Alliance healthcare	Sodium monohydrogenph...	12L05-N04	2
Alliance Healthcare	Aluminium potassium s...	13025404	1
Alliance Healthcare	Clotrimazole	AP1510003-01C	2
Alliance Healthcare	Lactose monohydrate	10195101	2
Alliance Healthcare	Macrogols 4000	12361512	1
Alliance Healthcare	Neomycin sulfate	1102043-03	6
Alliance Healthcare	Oxytetracycline hydro...	1205222-01	1
Alliance Healthcare	Oxytetracycline hydro...	1205028-02	15
Alliance Healthcare	Polymyxin B sulfate	1202039-02	5
Alliance Healthcare	Potassium sorbate	15382115	10
Alliance Healthcare	Progesterone, microni...	1206024-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Alliance Healthcare	Salicylic acid	12036215	5
Alliance Healthcare	Salicylic acid	12036216	6
Alliance Healthcare	Triclosan	13L-18-835	1
Alliance Healthcare/H...	Furosemide	12/13-N06	1
allianz	Testosterone propionate	15444704	1
Allianz Healthcare	Lidocaine	8565960	1
Allianz Healthcare	Pumpkin seed, peeled	16081210	1
Allianz Healthcare	Triamcinolone acetonide	16080005-01	1
almirall/gehe	Triamcinolone acetonide	17120607	8
alpha-cannabis Pharma...	Cannabidiol	ACPK911004	12
alpha-cannabis Pharma...	Cannabidiol	ACPK012405	1
Alte Apotheke	Erythromycin	16079206	5
Alte Apotheke	Lidocaine	16A28B06-321297	1
Alte Apotheke	Oak bark, cut	15314904	2
Alte Apotheke	Prednisolon dihydroge...	15K27B09-318738	1
Amos Vital	Ascorbic acid	L1120454	1
Amos Vital GmbH	Ascorbic acid	L1411003	1
Amos Vital GmbH/Allia...	Ascorbic acid	L1711004	1
an	Blessed thistle, cut	13389704	1
an	Camphor	12344909	1
an	Gentian root, cut	6311b130305	1
an	Hawthorn leaf with fl...	1971A130930	2
an	Passionflower herb, cut	13141212	1
an	Progesterone, microni...	1401002-01	3
an	Raspberry leaf, cut	13086708	1
an	Roselle	14018101	1
an	Sodium edetate	13142009	5
an	Triamcinolone acetonide	1403009-01	2
an	Verbena herb, cut	12383007	1
Andor	Metoprolol tartrate	APMTP0117MT	1
Andor	Testosterone propionate	1109024-03TP	1
Andreas Imhof	Menthol	8052015	1
Antares	Metoprolol tartrate	1014Q-07308	1
Anza	Benzoic acid	16001702	1
anzag	Prednisolone, microni...	14B18-B05-294966	13
anzag	Sorbitan stearate	14195	1
Anzag	Aluminium chloride he...	1204026-02	15
Anzag	Aluminium potassium s...	1218441	1
Anzag	Ambroxol hydrochloride	15G15-B01-316109	1
Anzag	Benzoic acid	16001708	1
Anzag	Betamethasone dipropi...	12114513	3
Anzag	Bilberry fruit, dried	2910N120921	2
Anzag	Bismuth gallath, basic	2012000059	5
Anzag	Blackberry leaf, cut	10136412	1
Anzag	Blessed thistle, cut	282281	1
Anzag	Calendula flower, cut	12221706	12
Anzag	Chloramphenicol	14346924	7
Anzag	Chlorhexidine diacetate	70558377	1
Anzag	Chlorhexidine diacetate	12177505	4
Anzag	Citric acid	10049712	1
Anzag	Citric acid	10049713	1
Anzag	Clobetasol propionate	15246001	7
Anzag	Clotrimazole	14C04-B30-2916618	1
Anzag	Clotrimazole	15G14-B06	2
Anzag	Clotrimazole	15395610	6
Anzag	Cocoa butter	12342004	3
Anzag	Dandelion herb with r...	274030	1
Anzag	Dexamethasone	17332404	9
Anzag	Dexamethasone	1233305	1
Anzag	Dithranol	120302502	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Anzag	Erythromycin	12152206	9
Anzag	Glucose	13329505	1
Anzag	Glucose monohydrate	17136805	13
Anzag	Glycerol monostearate...	12180502	1
Anzag	Hydrocortisone	11260305	1
Anzag	Hydrocortisone	1509229-01	1
Anzag	Hydroxyethylcellulose...	12249204	3
Anzag	Hydroxyethylcellulose...	1110024-03	23
Anzag	Isosorbide dinitrate ...	1510B01315217	1
Anzag	Ispaghula seed	11185818	1
Anzag	Lactose monohydrate	12139915	8
Anzag	Lavender flower	1488A12803	1
Anzag	Linseed	12137105	1
Anzag	Macrogols 4000	17068603	3
Anzag	Magnesium chloride he...	120802501	1
Anzag	Magnesium chloride he...	1304003-01	2
Anzag	Magnesium sulfate hep...	204035	1
Anzag	Maize style, cut	274920	1
Anzag	Meadowsweet herb, cut	10381604	5
Anzag	Menthol	16L20-T04-044769	10
Anzag	Menthol	12/19-N10	2
Anzag	Metronidazole	11373512	1
Anzag	Miconazole nitrate	12374301	8
Anzag	Miconazole nitrate	15A13-B07-308432	3
Anzag	Minoxidil	13E13-N10	12
Anzag	Mometasone furoate, m...	12J29N03	1
Anzag	Nettle leaf, cut	2381-A130122	1
Anzag	Nystatin	2902041-02	1
Anzag	Oxytetracycline hydro...	1506037	4
Anzag	Passionflower herb, cut	1207211669	1
Anzag	Prednisolone acetate	13A07N07	1
Anzag	Prednisolone acetate	16246812	5
Anzag	Procaine hydrochloride	60998058	1
Anzag	Progesterone, microni...	15C31-B06-312963	5
Anzag	Progesterone, microni...	2515A-03235	1
Anzag	Progesterone, microni...	14k11-b04-300921	1
Anzag	Progesterone, microni...	13K28-B01-293478	1
Anzag	Propyl-4-hydroxybenzo...	11373512	5
Anzag	Raspberry leaf, cut	2291A120806	2
Anzag	Raspberry leaf, cut	12103107	11
Anzag	Sodium citrate	1512705	1
Anzag	Sodium edetate	13K19-B02-290051	9
Anzag	Sweet woodruff herb, ...	12103504	2
Anzag	Tetracaine hydrochlor...	17256705	2
Anzag	Tretinoin	15078305	3
Anzag	Trichloroacetic acid	12135204	2
Anzag	Triclosan	5846205	4
Anzag	Trometamol	12113512	13
Anzag	Wax, yellow	1462	1
Anzag	Xylometazoline hydroc...	12027911	1
ANZAG	Progesterone, microni...	12195907	6
ANZAG	Walnut leaf, cut	12124603	11
Anzag 21.10.14 EK 2...	Metronidazole	13095007	6
Anzag / Noweda	Benzocaine	1321351	1
Anzag 05.03.13 ek: 9,...	Wheat starch	104A120621	3
Anzag 05.12.12	Agrimony, cut	11374002	3
Anzag 06.04.2016	Aniseed	15398601	1
Anzag 10.05.13 ek:3,16EUR	Dandelion herb with r...	6611B130321	1
Anzag 12.03.13 ek: 4,...	St. John's wort, cut	11185608	2
Anzag 12.11.12 ek:3,85	Yarrow herb, cut	11280513	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Anzag 14.11.12 ek:40,14	Minoxidil	120303501	1
Anzag 14.12.12 ek:10,01	Hydrocortisone	1202032-01	5
Anzag 17.04.13 ek: 2,...	Tall golden rod, cut	4521A110816	2
anzag 18.12.14 25,19	Methacholine chloride	14G30-B03-296209	1
anzag 20.11.2013	Triclosan	12359114	3
Anzag 22.03.13 ek: 5,...	Aloes, cape, powdered	12232102	1
ANZAG 23.07.14	Ketoconazole	13h09-n07	1
Anzag 23.11.12 ek:26,49	Erythromycin	1205035-01	15
Anzag 25.10.12 ek:3.16EUR	Dandelion herb with r...	6611B110131	1
Anzag 26.10.12 2.36EUR	Oak bark, cut	10348616	6
Anzag 27.12.12 ek:2,79	Blackberry leaf, cut	10136048	1
Anzag 29.04.13 ek: 4,...	Magnesium sulfate hep...	12138518	3
Anzag ek: 15,60 20.07...	Iceland moss, cut	11101605	1
Anzag ek: 2,59EUR 19.09...	Caraway fruit	11003216	3
Anzag ek: 3,06 11.09...	Sage leaf, cut	2301A120403	38
Anzag ek: 4,58EUR 25.09...	Walnut leaf, cut	2111A120123	8
Anzag ek: 4,88EUR 26.09...	Male speedwell wort, ...	4661C120321	10
Anzag ek: 6,39 11.09...	Alchemilla, cut	11182806	6
Anzag ek:2,50EUR 15.09.12	Horse chestnut leaf, ...	2551B111013	7
Anzag ek:2,64EUR 16.07.12	Mistletoe, cut	10273417	1
Anzag ek:2,98EUR 25.08.12	Chloramphenicol	11370711	13
Anzag ek:3,00EUR 11.09.12	Raspberry leaf, cut	12103105	18
Anzag ek:3,13EUR 15.09.12	Melilot herb, cut	11307605	4
Anzag ek:7,89EUR	Dequalinium chloride	11166204	3
Anzag ek:8,77EUR 30.08.12	Manna, cut	11308904	2
Anzag EK4,31 17.08.12	Nettle leaf, cut	2381A110913	19
Anzag EK7,91 08.08.12	Gentian root, cut	10263407	1
Anzag/ Hedinger	Sodium hydrogen carbo...	206027	1
Anzag/Anzag 11.09.18	Riboflavin	18B22-B01-347142	1
Anzag/Anzag 11.09.18	Salicylic acid	18F05-B02-350586	3
Anzag/Anzag 26.07.18	Calcium pantothenate	1803017-01	1
Anzag/euRho	Tetracycline hydrochl...	1110029-01Nachprüfun	1
Anzag/Fragon	Lactose monohydrate	1005SLP	1
Anzag/Henry Lamotte	Salicylic acid	1334018	1
Anzag/Henry Lamotte	Tretinoin	15078303	1
Anzag/Henry Lamotte	Urea	13-18-B03-303581	1
Anzag: 07.09.12 ek:4.01	Chlortetracycline hyd...	12C16-N01	1
Anzag;18.03.15;2,76EUR	Allantoin	14247602	8
Apo Ident/ Noweda	Gentamicin sulfate	1711025	1
Apocare	Cannabidiol	ACPK312001	1
ApoFit	Blue bottle flower wi...	278776	1
ApoFit	Fenugreek	222392	1
ApoFit/Noweda	Walnut leaf, cut	271137	1
apoident	Triamcinolone acetone...	1508008-01	1
apo-ident.	Minoxidil	1607007-07	1
apomix	Tartaric acid		2
apomix	Tartaric acid	1111	1
Apomix	Prednisolone, microni...	16H04-B03-332930	6
Apomix	Trichloroacetic acid	12G04-M11	1
Apomix / 11	Carbachol	12330703	1
Apomix / AHD	Testosterone propionate	40/14	1
Apomix/Kehr	Oxytetracycline hydro...	16077014	1
Apomix/Noweda	Octenidine dihydrochl...	13166813	1
Apomix/Phönix	Triamcinolone acetone...	13k15-b31-289029	2
APOrtha Deutschland G...	Magnesium citrate	17000809/0	1
Apotheke Bühlau	Almonds, bitter, whole	267467	2
Apotheke Bühlau	Horse chestnut leaf, ...	100416	1
Apotheke im Forum	Glucose monohydrate	17010002-01	1
Apotheke/Defektur Mis...	Glucose	39/15(02.11.2015)	1
Apotheke-SiO2-Mischung	Glucose	18/15	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Arktis	Calcium ascorbate	7339/4617	1
Arktis	Calcium ascorbate	733P/1518	1
Arktis	Magnesium citrate	894P/2418	1
Ashland Ind/Gehe	Hydroxyethylcellulose...	140017-01	1
Ashland Industries/An...	Hydroxyethylcellulose...	L1891	1
Ashland/Sanacorp	Neomycin sulfate	1303033-02	15
Audor Pharma	Alfatradiol	AP5100461A	1
Audor Pharma	Alfatradiol	33-16-I-62P2985A	1
Audor Pharma	Alfatradiol	AP33-16-I-42P2945A	3
Audor Pharma	Alfatradiol	APXPC0071A	3
Audor Pharma	Alfatradiol	AP33-18-I-68P3620A	1
Audor Pharma	Alfatradiol	AP5100461A	1
Audor Pharma	Alfatradiol	APXP103DAA	1
Audor Pharma	Alfatradiol	AP33-16-I-62P2985A	3
Audor Pharma	Alfatradiol	AP1404001-01A	1
Audor Pharma	Alfatradiol	APXP103DAA	1
Audor Pharma	Alfatradiol	AP5104925A	2
Audor Pharma	Alfatradiol	AP5104925A	4
Audor Pharma	Alfatradiol	AP21-16-II-14P4511A	4
Audor Pharma	Alfatradiol	AP1504015-01A	2
Audor Pharma	Alfatradiol	AP514925A	1
Audor Pharma	Alfatradiol	AP3316I62P2985A	1
Audor Pharma	Benzydamine hydrochlo...	300111702007B	1
Audor Pharma	Benzydamine hydrochlo...	AP300111702007B	3
Audor Pharma	Benzydamine hydrochlo...	300111702007b	1
Audor Pharma	Benzydamine hydrochlo...	AP300111801004B	1
Audor Pharma	Citric acid	PA201	1
Audor Pharma	Clobetasol propionate	APL035E12C	1
Audor Pharma	Clobetasol propionate	APCL0102/0518C	1
Audor Pharma	Clobetasol propionate	sieheBemerkung	1
Audor Pharma	Clobetasol propionate	APB0054/112010C	1
Audor Pharma	Clobetasol propionate	APL/035/E-12C	1
Audor Pharma	Clobetasol propionate	APCLP/M/036/16C	1
Audor Pharma	Clotrimazole	APCLT/1703/0016CC	2
Audor Pharma	Clotrimazole	APCLT/1703/0016C	1
Audor Pharma	Clotrimazole	1703/0016CC	1
Audor Pharma	Clotrimazole	AP1306036-01C	1
Audor Pharma	Colistin sulfate	APKB-CS-1603034C	2
Audor Pharma	Colistin sulfate	APKB-CS-1610012C	1
Audor Pharma	Colistin sulfate	APKB-CS-1601026C	1
Audor Pharma	Cyproterone acetate	APCPA140201CP	1
Audor Pharma	Cyproterone acetate	cpa140201cp	1
Audor Pharma	Cyproterone acetate	X03-181002-m0C	1
Audor Pharma	Dithranol	APOP-DNL/11/17/001	1
Audor Pharma	Dithranol	APDT/375/2012D	11
Audor Pharma	Dithranol	APDT3902013D	1
Audor Pharma	Dithranol	APOP-DNL11170011D	1
Audor Pharma	Dithranol	OP-DNL/11/17/001-1D	1
Audor Pharma	Dithranol	DT/351/2010D	1
Audor Pharma	Dithranol	DT/375/2012D	3
Audor Pharma	Dithranol	APOP-DNL/11/17/001D	1
Audor Pharma	Dithranol	AP0-DNL/11/17/001-1D	1
Audor Pharma	Dithranol	APOP-DNL/11/17/001-	2
Audor Pharma	Dithranol	APDT/365/2011D	1
Audor Pharma	Dithranol	28 (APOP-DNL/11/17/0	1
Audor Pharma	Dithranol	APDT/365/2011	2
Audor Pharma	Dithranol	A00P-DNL/11/17/001-	1
Audor Pharma	Dithranol	APDT/351/2010D	1
Audor Pharma	Dithranol	APOP-DNL/11/17/001-1D	1
Audor Pharma	Dithranol	AP-OP-DNL/11/17/001-	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Audor Pharma	Emulsifying cetostear...	16238602	1
Audor Pharma	Erythromycin	EB022I17E	1
Audor Pharma	Estriol	AP207800.41.02.0217E	1
Audor Pharma	Estriol	AP20780041020218E	1
Audor Pharma	Estriol	APL00037293	1
Audor Pharma	Estriol	AP207800.41.03.0616	1
Audor Pharma	Estriol	APL000036940E	1
Audor Pharma	Estriol	AP207800.41.03.0616E	2
Audor Pharma	Estriol	APL00036940E	2
Audor Pharma	Estriol	APL00037293E	1
Audor Pharma	Glucose	16157402	6
Audor Pharma	Glycerol monostearate...	13232604	4
Audor Pharma	Hyaluronic acid sodiu...	AP30101305H	2
Audor Pharma	Hyaluronic acid sodiu...	AP30101505H	1
Audor Pharma	Hyaluronic acid sodiu...	AP30101404H	1
Audor Pharma	Hyaluronic acid sodiu...	AP30101404H	1
Audor Pharma	Melatonin	ALC/MLT/170901	1
Audor Pharma	Methoxsalen	AP1011024-01M	1
Audor Pharma	Methoxsalen	AP15244213M	6
Audor Pharma	Methoxsalen	1203009M	1
Audor Pharma	Methoxsalen	APHCFMPC16011M	7
Audor Pharma	Methoxsalen	15244213	1
Audor Pharma	Metoprolol tartrate	APMTPOC117MT	1
Audor Pharma	Metoprolol tartrate	APMTPO2218MT	1
Audor Pharma	Metoprolol tartrate	MTP02218	1
Audor Pharma	Metoprolol tartrate	APMTPO0117MT	7
Audor Pharma	Metoprolol tartrate	MTP00117MT	1
Audor Pharma	Minoxidil	AP-016/15M	1
Audor Pharma	Minoxidil	K007.1/16/04/065-1	1
Audor Pharma	Minoxidil	AP016/15M	7
Audor Pharma	Minoxidil	AP007/15M	3
Audor Pharma	Minoxidil	APK007	1
Audor Pharma	Minoxidil	8/14M	1
Audor Pharma	Minoxidil	APK007.1/16/04/065	2
Audor Pharma	Minoxidil	K007.1/16/04/065M	1
Audor Pharma	Minoxidil	AP008/14M	2
Audor Pharma	Minoxidil	APK007.1/16/04/065-1	2
Audor Pharma	Minoxidil	AP007/15M	1
Audor Pharma	Minoxidil	APK007.1/16/04/065M	6
Audor Pharma	Minoxidil	AP012/16M/5270	1
Audor Pharma	Minoxidil	AP020/11M	3
Audor Pharma	Minoxidil	414E-02844	1
Audor Pharma	Minoxidil	AP006/17M	3
Audor Pharma	Minoxidil	AP012/16M	4
Audor Pharma	Minoxidil	APK007.1/17/12/191M	6
Audor Pharma	Minoxidil	AP010/12M	5
Audor Pharma	Minoxidil	K007.1/17/12/191M	2
Audor Pharma	Minoxidil	K007.1/16/04/065-1M	1
Audor Pharma	Minoxidil	APK007.1/17/09/139M	5
Audor Pharma	Minoxidil	APK007.1/17/12/191	2
Audor Pharma	Minoxidil	AP001/16M	7
Audor Pharma	Minoxidil	AP088/14M	1
Audor Pharma	Minoxidil	AP008/14M	4
Audor Pharma	Minoxidil	APK007.1/16/04/065-	3
Audor Pharma	Mometasone furoate, m...	APMOM006/0177M	1
Audor Pharma	Mometasone furoate, m...	APMOM301/0118M	2
Audor Pharma	Mometasone furoate, m...	APMOM006/0117M	4
Audor Pharma	Mometasone furoate, m...	APMOM006/0517M	5
Audor Pharma	Naloxone hydrochlorid...	AP13501078N	1
Audor Pharma	Naphazoline nitrate	12I21-N05	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Audor Pharma	Oxybutynin hydrochlor...	AP0H/043/2017-2018	1
Audor Pharma	Oxybutynin hydrochlor...	AP0h/012/2017-2018-	1
Audor Pharma	Oxybutynin hydrochlor...	AP0H0122017-2018-10B	1
Audor Pharma	Oxybutynin hydrochlor...	AP0H/012/2017-2018-	1
Audor Pharma	Oxybutynin hydrochlor...	AP0H0122017201810B	2
Audor Pharma	Prasterone	AP11040002P	2
Audor Pharma	Prasterone	10050051/1P	1
Audor Pharma	Prasterone	APZDHA16003-M	1
Audor Pharma	Prasterone	APZDHAC17003-MD	1
Audor Pharma	Pregnenolone	APC07-150801(M)P	3
Audor Pharma	Pregnenolone	APC07-180201(M)P	3
Audor Pharma	Pregnenolone	AP10L10-M01P	1
Audor Pharma	Pregnenolone	AP1109007-01P	1
Audor Pharma	Pregnenolone	APC07-180201	2
Audor Pharma	Pregnenolone	APC07-160801(M)P	2
Audor Pharma	Progesterone, microni...	APC01-160201P	2
Audor Pharma	Progesterone, microni...	APC01-160201	1
Audor Pharma	Progesterone, microni...	APC01-181003P	11
Audor Pharma	Progesterone, microni...	APC01-171101P	31
Audor Pharma	Progesterone, microni...	APC01-170828P	1
Audor Pharma	Progesterone, microni...	APC01-151114P	3
Audor Pharma	Progesterone, microni...	APC01-180610P	10
Audor Pharma	Progesterone, microni...	AP0C01-1700828P	2
Audor Pharma	Progesterone, microni...	C01-180128	10
Audor Pharma	Progesterone, microni...	C01-120708P	1
Audor Pharma	Progesterone, microni...	APC01-150311P	1
Audor Pharma	Progesterone, microni...	APC01-120613P	1
Audor Pharma	Progesterone, microni...	APC01-120708P	3
Audor Pharma	Progesterone, microni...	APC01-141212P	1
Audor Pharma	Progesterone, microni...	APC01-150708P	4
Audor Pharma	Sodium citrate	171000601	1
Audor Pharma	Sucrose	16206102	13
Audor Pharma	Testosterone	apztstg18004t	2
Audor Pharma	Testosterone	APZTST16001-M1T	1
Audor Pharma	Testosterone	APZTSTG18004T	2
Audor Pharma	Testosterone propionate	AP12J10-N08TP	2
Audor Pharma	Testosterone propionate	AP1109024-03	1
Audor Pharma	Testosterone propionate	AP13K25-B05-287643T	1
Audor Pharma	Testosterone propionate	AP12E04N01TP	2
Audor Pharma	Testosterone propionate	APZTSP1503-MTP	1
Audor Pharma	Testosterone propionate	APZTSP15003-MTP	11
Audor Pharma	Testosterone propionate	AP1109024-03	1
Audor Pharma	Testosterone propionate	AP13K25-B05-287643TP	1
Audor Pharma	Testosterone propionate	AP1109024-03TP	4
Audor Pharma	Tranexamic acid	1508016MT	1
Audor Pharma	Tranexamic acid	APX1508016MT	1
Audor Pharma	Tranexamic acid	X1712011M	1
Audor Pharma	Tranexamic acid	APX17112011MT	1
Audor Pharma	Tranexamic acid	X1611013M	2
Audor Pharma	Tranexamic acid	APX1611013MT	2
Audor Pharma	Tranexamic acid	APX1712011MT	6
Audor Pharma	Triamcinolone acetonide	APU0544/316270T	2
Audor Pharma	Triamcinolone acetonide	15240T	1
Audor Pharma	Triamcinolone acetonide	APU0544/315240T	2
Audor Pharma	Triamcinolone acetonide	APU0544/3162705	1
Audor Pharma	Triamcinolone acetonide	APU0544/31524T	1
Audor Pharma	Triamcinolone acetonide	APU0544/317171T	1
Audor Pharma	Triclosan	APK001.1/17/03/139T	2
Audor Pharma	Ubidecarenone	APSCFQ150614C	1
Audor Pharma	Ubidecarenone	APSZFQ150614C	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Audor Pharma	Urea	PU145	1
Audor Pharma	Uridine	APIB1400280AU	1
Audor Pharma	Uridine	AP721060040/5U	3
Audor Pharma	Uridine	AP30514U	4
Audor Pharma	Uridine	APIB1600182.11AU	1
Audor Pharma	Uridine	APIB1600182AU	6
Audor Pharma	Uridine	AP30514IB1500001AXU	4
Audor Pharma	Uridine	APIB1600169AU	22
Audor Pharma	Uridine	APIB1600169	1
Aug. Hedinger	Sodium citrate	8765456	1
Aurica	Aniseed	521302895	1
Aurica	Birch leaf, cut	521205644	1
Aurica	Common wormwood herb,...	521305386	1
Aurica	Dandelion herb with r...	521200197	1
Aurica	Red clover flower	521206331	1
Aurica	Red clover flower	521203442	1
Aurica	Rose hips with peel, ...	521201801	1
Aurica	Sage leaf, cut	521204613	1
Aurica	St. John's wort, cut	521200778	1
Aurica	St. John's wort, cut	521202559	1
Aurica	Tall golden rod, cut	521201267	1
Aurica	Tall golden rod, cut	224133	1
Aurica	Thyme, rubbed	521207003	1
Aurica	Valerian root, cut	221838	2
Aurica	Verbena herb, cut	521305668	1
Aurica / 11	Equisetum stem, cut	521302898	1
Aurica / 11	Juniper	521301390	1
Aurica / 11	Primrose flower	521301234	1
Aurica Naturheilmitte...	Aniseed, powdered	521405796	1
Aurica/Anzag	Red clover flower	521303918	1
Aurica/Gehe	Primrose flower	521203308	1
aurica/phoenix	Liquorice root, natur...	521206925	1
Aurica/Phönix	Acorus root, peeled	150649	1
Aurica/Phönix	Angelica root	520800638	1
Aurica/Phönix	Bean shells, cut	5052301	1
Aurica/Phönix	Centaury, cut	520702152	1
Aurica/Phönix	Juniper	521006043	1
Aurica/Phönix	Maize style, cut	520901387	1
Aurica/Phönix	Male speedwell wort, ...	520708362	1
Aurica/Phönix	Mistletoe, cut	6110801	1
Axyntis Group	Clioquinol	1301014-03	1
Bahnhof-Apotheke	Acorus root, peeled	17-121107	1
Bahnhof-Apotheke	Alchemilla, cut	17-121107	1
Bahnhof-Apotheke	Aniseed	7-120913	1
Bahnhof-Apotheke	Birch leaf, cut	18-120702	1
Bahnhof-Apotheke	Bitter fennel	11-121029	1
Bahnhof-Apotheke	Blackberry leaf, cut	18-120702	1
Bahnhof-Apotheke	Calendula flower, cut	17-121106	1
Bahnhof-Apotheke	Caraway fruit	18-120709	1
Bahnhof-Apotheke	Elder flower, rubbed	17-121112	1
Bahnhof-Apotheke	Hawthorn leaf with fl...	7-121011	1
Bahnhof-Apotheke	Hay flower	13-121105	1
Bahnhof-Apotheke	Lactose monohydrate	120614	1
Bahnhof-Apotheke	Lactose monohydrate	120808	1
Bahnhof-Apotheke	Lactose monohydrate	121119	2
Bahnhof-Apotheke	Lactose monohydrate	130115	1
Bahnhof-Apotheke	Lactose monohydrate	130729	1
Bahnhof-Apotheke	Lactose monohydrate	130911	1
Bahnhof-Apotheke	Lactose monohydrate	130925	1
Bahnhof-Apotheke	Lactose monohydrate	131023	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bahnhof-Apotheke	Lactose monohydrate	131112	1
Bahnhof-Apotheke	Lactose monohydrate	131127	2
Bahnhof-Apotheke	Lactose monohydrate	140630	1
Bahnhof-Apotheke	Lactose monohydrate	131218	2
Bahnhof-Apotheke	Lactose monohydrate	140113	3
Bahnhof-Apotheke	Lactose monohydrate	140220	1
Bahnhof-Apotheke	Lactose monohydrate	140206	2
Bahnhof-Apotheke	Lactose monohydrate	140310	2
Bahnhof-Apotheke	Lime flower, cut	7-121011	1
Bahnhof-Apotheke	Liquorice root, peeled	17-121112	1
Bahnhof-Apotheke	Lycopodium herb, cut	17-121112	1
Bahnhof-Apotheke	Maize style, cut	20-121018	1
Bahnhof-Apotheke	Male speedwell wort, ...	4-121002	1
Bahnhof-Apotheke	Matricaria flower	11-121024	1
Bahnhof-Apotheke	Melissa leaf, cut	17-120917	1
Bahnhof-Apotheke	Nettle leaf, cut	7-121109	1
Bahnhof-Apotheke	Oak bark, cut	17-121106	1
Bahnhof-Apotheke	Psyllium seed	3-121005	1
Bahnhof-Apotheke	Raspberry leaf, cut	17-120917	1
Bahnhof-Apotheke	Ribwort plantain, cut	7-121011	1
Bahnhof-Apotheke	Rosemary leaf, cut	7-120915	1
Bahnhof-Apotheke	Sage leaf, cut	3-121018	1
Bahnhof-Apotheke	Shepherd's purse, cut	4-120924	1
Bahnhof-Apotheke	St. John's wort, cut	4-121018	1
Bahnhof-Apotheke	Thyme, rubbed	17-121112	1
Bahnhof-Apotheke	Willowherb, small flo. ...	17-121107	1
Bahnhof-Apotheke	Wormwood herb, cut	17-121107	1
Bahnhof-Apotheke	Yarrow herb, cut	7-121011	1
Bahnhofsapotheke	Lactose monohydrate	140401	1
BASF Grenzach GmbH/Ph...	Triclosan	13A11-No3	1
BASF Grenzach/Noweda	Triclosan	5900346	1
BASF/Phönix	Tretinoin	1304028-02	8
Bayer Pharma AG	Estradiol hemihydrate	1411024-01	4
Bayer-Schering	Estradiol benzoate	1007037-03	1
Beiersdorf	Triclosan	14F23-B02-295764	3
Beiersdorf AG	Tartaric acid	11G19-N08	1
Beiersdorf/Phönix	Glucose monohydrate	L6209	1
Berkel AHK	Agnus castus fruit	14194401	4
Berkel AHK	Aniseed	11374103	1
Berkel AHK	Aniseed	14194401	1
Berkel AHK	Calcium gluconate	14085407	2
Berkel AHK	Oak bark, cut	13360310	1
Berkel AHK	Prednisolone, microni...	13343115	2
Berkel AHK	Valerian root, powdered	13287406	1
Berkel AHK	Valerian root, powdered	14179403	2
Biomed	Erythromycin	1508035-01	1
Bionorica	Ascorbyl palmitate	85532	5
Bionorica	Ascorbyl palmitate	74736	2
Bionorica	Ascorbyl palmitate	74962	1
Bionorica /Phönix	Miconazole nitrate	12C13B03-279593	1
Bionorica Ethics	Ascorbyl palmitate	82275	4
Bionorica/Noweda	Ascorbyl palmitate	77773	2
Bios	Lu Hui	424263	1
bohn	Savory, rubbed	12177208	1
Bomastus / Alliance H...	Miconazole nitrate	16E04-B03-327367	1
Bomastus-Werke/Phönix	Rosemary leaf, cut	298424	2
Bomba/AHD	Equisetum stem, cut	280449	3
Bomba/AHD	Senna leaf, cut	280008	1
Bomba/AHD	Tall golden rod, cut	279526	1
Bomba/AHD	Willow bark, cut	280362	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Hop strobile	274724	1
Bombastus	Acorus root, natural,...	276260	2
Bombastus	Acorus root, natural,...	274159	1
Bombastus	Acorus root, natural,...	224114	1
Bombastus	Acorus root, natural,...	281692	1
Bombastus	Acorus root, natural,...	277958	1
Bombastus	Acorus root, peeled	274159	3
Bombastus	Acorus root, peeled	267365	1
Bombastus	Alchemilla, cut	276173	3
Bombastus	Alchemilla, cut	14066007	1
Bombastus	Alchemilla, cut	231013B	1
Bombastus	Alchemilla, cut	521203995	1
Bombastus	Alchemilla, cut	282032	1
Bombastus	Alchemilla, cut	278889	1
Bombastus	Alchemilla, cut	266404	1
Bombastus	Alchemilla, cut	279070	4
Bombastus	Alchemilla, cut	272021	5
Bombastus	Almonds, bitter, whole	283369	1
Bombastus	Almonds, bitter, whole	303218	2
Bombastus	Almonds, bitter, whole	304961	1
Bombastus	Almonds, bitter, whole	1700352	2
Bombastus	Almonds, bitter, whole	267529	1
Bombastus	Almonds, bitter, whole	284978	1
Bombastus	Almonds, bitter, whole	275001	2
Bombastus	Almonds, bitter, whole	275601	1
Bombastus	Almonds, bitter, whole	297921	1
Bombastus	Almonds, bitter, whole	275600	1
Bombastus	Almonds, bitter, whole	305322	2
Bombastus	Almonds, bitter, whole	306048	1
Bombastus	alpha-tocopheryl acet...	1507001-02	1
Bombastus	Angelica root	305460	1
Bombastus	Angelica root	272895	1
Bombastus	Angelica root	296334/1121115	1
Bombastus	Angelica root	301768	1
Bombastus	Angelica root	291848	1
Bombastus	Angelica root	274128	2
Bombastus	Angelica root	295085	2
Bombastus	Angelica root	269775	2
Bombastus	Angelica root	272954	1
Bombastus	Angelica root	275879	1
Bombastus	Angelica root	275878	1
Bombastus	Angelica root	270117	1
Bombastus	Angelica root	278757	4
Bombastus	Angelica root	281280	7
Bombastus	Angelica root	278294	1
Bombastus	Angelica root	307041	3
Bombastus	Angelica root	290309	1
Bombastus	Aniseed	305177	1
Bombastus	Aniseed	285778	2
Bombastus	Aniseed	282917	2
Bombastus	Aniseed	285777	1
Bombastus	Aniseed	289846	2
Bombastus	Artichoke leaf, cut	8202213593	1
Bombastus	Ascorbyl palmitate	525526	1
Bombastus	Bean shells, cut	282475	1
Bombastus	Bean shells, cut	278686	1
Bombastus	Bean shells, cut	278696	2
Bombastus	Bean shells, cut	240503	1
Bombastus	Bean shells, cut	270517	1
Bombastus	Bean shells, cut	272781	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Bearberry leaf, cut	275491	2
Bombastus	Bearberry leaf, cut	267926	1
Bombastus	Bearberry leaf, cut	272297	8
Bombastus	Betamethasone dipropi...	12B27N08	1
Bombastus	Betamethasone valerate	1510037-01	23
Bombastus	Betony herb, cut	12007205	2
Bombastus	Bilberry leaf, cut	272176	2
Bombastus	Birch leaf, cut	272385	1
Bombastus	Birch leaf, cut	287280	1
Bombastus	Birch leaf, cut	280677	1
Bombastus	Birch leaf, cut	270743	1
Bombastus	Birch leaf, cut	267077	1
Bombastus	Bitter fennel	285609	1
Bombastus	Bitter fennel	289870	2
Bombastus	Bitter fennel	289869	2
Bombastus	Bitter fennel	275922	4
Bombastus	Bitter fennel	299386	3
Bombastus	Bitter fennel	296877	3
Bombastus	Bitter fennel	268130	2
Bombastus	Bitter fennel	289870	1
Bombastus	Bitter fennel	299385	2
Bombastus	Bitter fennel	293627	1
Bombastus	Bitter fennel	280891	2
Bombastus	Bitter orange peel, cut	278233	4
Bombastus	Bitter orange peel, cut	274240	3
Bombastus	Blackberry leaf, cut	267122	1
Bombastus	Blackberry leaf, cut	273668	2
Bombastus	Blessed thistle, cut	271946	1
Bombastus	Blue bottle flower wi...	275759	2
Bombastus	Blue bottle flower wi...	284030	1
Bombastus	Calendula flower	280306	1
Bombastus	Calendula flower	268626	1
Bombastus	Calendula flower, cut	275813	1
Bombastus	Calendula flower, cut	234997	1
Bombastus	Calendula flower, cut	253620	1
Bombastus	Calendula flower, cut	270046	1
Bombastus	Calendula flower, cut	207294	1
Bombastus	Calendula flower, cut	279632	1
Bombastus	Calendula flower, cut	275128	1
Bombastus	Caraway fruit	275849	1
Bombastus	Caraway fruit	270124	1
Bombastus	Caraway fruit	282269	1
Bombastus	Carmellose sodium (40...	13104813	1
Bombastus	Centaury, cut	277598	3
Bombastus	Cetostearyl alcohol	295148	1
Bombastus	Cetostearyl alcohol	300132	2
Bombastus	Citric acid	272757	2
Bombastus	Citric acid monohydrate	294630	1
Bombastus	Citric acid monohydrate	309650	1
Bombastus	Citric acid monohydrate	309344	2
Bombastus	Citric acid monohydrate	70620101	1
Bombastus	Citric acid monohydrate	309343	1
Bombastus	Citric acid monohydrate	259437	1
Bombastus	Clobetasol propionate	13G23-N0	1
Bombastus	Clotrimazole	Validierungsmessung	1
Bombastus	Coltsfoot leaf, cut	277848	2
Bombastus	Coltsfoot leaf, cut	274289	5
Bombastus	Coltsfoot leaf, cut	276303	5
Bombastus	Coltsfoot leaf, cut	281100	4
Bombastus	Coltsfoot leaf, cut	272019	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Coltsfoot leaf, cut	269501	3
Bombastus	Coltsfoot leaf, cut	275489	3
Bombastus	Coltsfoot leaf, cut	265856	1
Bombastus	Coltsfoot leaf, cut	278758	7
Bombastus	Coltsfoot leaf, cut	285263	3
Bombastus	Coltsfoot leaf, cut	277131	1
Bombastus	Coltsfoot leaf, cut	232339	1
Bombastus	Comfrey herb, cut	281531	1
Bombastus	Coriander	263733	1
Bombastus	Couch grass rhizome, ...	262175	1
Bombastus	Couch grass rhizome, ...	276777	3
Bombastus	Couch grass rhizome, ...	273387	3
Bombastus	Couch grass rhizome, ...	270500	2
Bombastus	Couch grass rhizome, ...	278926	1
Bombastus	Dandelion herb with r...	282916	7
Bombastus	Dandelion herb with r...	274030	4
Bombastus	Dandelion herb with r...	280363	5
Bombastus	Dandelion herb with r...	290125	1
Bombastus	Dandelion herb with r...	285264	1
Bombastus	Dandelion herb with r...	277597	6
Bombastus	Dandelion herb with r...	271272	5
Bombastus	Dandelion herb with r...	258809	1
Bombastus	Dandelion herb with r...	280363-3385310110	1
Bombastus	Dandelion herb with r...	275126	1
Bombastus	Dexamethasone	150546629	1
Bombastus	Dill herb, cut	50836505	1
Bombastus	Elder flower, rubbed	269873	2
Bombastus	Elder flower, rubbed	275393	4
Bombastus	Elder flower, rubbed	265631	1
Bombastus	Elder flower, rubbed	279216	3
Bombastus	Emulsifying cetostear...	308514	1
Bombastus	Emulsifying cetostear...	265360	1
Bombastus	Emulsifying cetostear...	295148	1
Bombastus	Emulsifying cetostear...	303935	4
Bombastus	Emulsifying cetostear...	308514	2
Bombastus	Emulsifying cetostear...	288691	3
Bombastus	Emulsifying cetostear...	300132	4
Bombastus	Emulsifying cetostear...	18003404	5
Bombastus	Emulsifying cetostear...	278418	1
Bombastus	Emulsifying cetostear...	18003437	3
Bombastus	Equisetum stem, cut	6011401	1
Bombastus	Equisetum stem, cut	268321	1
Bombastus	Equisetum stem, cut	275891	3
Bombastus	Equisetum stem, cut	265924	1
Bombastus	Equisetum stem, cut	274714	1
Bombastus	Equisetum stem, cut	280647	1
Bombastus	Equisetum stem, cut	272552	2
Bombastus	Equisetum stem, cut	278107	1
Bombastus	Equisetum stem, cut	273963	1
Bombastus	Equisetum stem, cut	264101	1
Bombastus	Equisetum stem, cut	264100	1
Bombastus	Equisetum stem, cut	18111304	1
Bombastus	Estriol	5128	1
Bombastus	Euphrasy herb, cut	274553	1
Bombastus	Euphrasy herb, cut	270726	2
Bombastus	Euphrasy herb, cut	278108	2
Bombastus	Euphrasy herb, cut	277846	1
Bombastus	Frangula bark, cut	275998	2
Bombastus	Frangula bark, cut	268407	1
Bombastus	Frangula bark, cut	273926	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Frangula bark, cut	270733	1
Bombastus	Frangula bark, cut	281284	2
Bombastus	Frangula bark, cut	278599	3
Bombastus	Frangula bark, cut	238812	1
Bombastus	Gentamicin sulfate	700175	1
Bombastus	Gentamicin sulfate	13024-N05	1
Bombastus	Ginseng root	13150903	1
Bombastus	Glucose	160505202	1
Bombastus	Griseofulvin	15L03-HH150640	1
Bombastus	Hard fat in pastilles	302709	1
Bombastus	Hard fat in pastilles	293848	6
Bombastus	Hard fat in pastilles	306363	4
Bombastus	Hard fat in pastilles	308559	1
Bombastus	Hard fat in pastilles	18000346	4
Bombastus	Hard fat in pastilles	287287	1
Bombastus	Hard fat in pastilles	273759	1
Bombastus	Hard fat W35	302709	2
Bombastus	Hawthorn leaf with fl. ...	272881	1
Bombastus	Hawthorn leaf with fl. ...	267551	1
Bombastus	Hawthorn leaf with fl. ...	264044	1
Bombastus	Hawthorn leaf with fl. ...	275310	2
Bombastus	Hawthorn leaf with fl. ...	273287	1
Bombastus	Hawthorn leaf with fl. ...	278756	1
Bombastus	Hawthorn leaf with fl. ...	267000	1
Bombastus	Heather, cut	278314	1
Bombastus	Hop strobile	276976	1
Bombastus	Hop strobile	271870	1
Bombastus	Hop strobile	208502	1
Bombastus	Hydrocortisone	14I05-B02-300227	1
Bombastus	Hydrocortisone acetate	2712M01718	1
Bombastus	Hydrocortisone acetate	123Validierung	1
Bombastus	Iceland moss, cut	13357802	1
Bombastus	Lactose monohydrate	305227	2
Bombastus	Lactose monohydrate	306147	3
Bombastus	Lactose monohydrate	305336/0	1
Bombastus	Lactose monohydrate	292983	1
Bombastus	Lactose monohydrate	286354	1
Bombastus	Lactose monohydrate	286710	1
Bombastus	Lactose monohydrate	292995	1
Bombastus	Lactose monohydrate	283992	2
Bombastus	Lactose monohydrate	282648	1
Bombastus	Lactose monohydrate	308118/0	1
Bombastus	Lactose monohydrate	307702/0	1
Bombastus	Lactose monohydrate	307976/0	1
Bombastus	Lactose monohydrate	306651/0	1
Bombastus	Lactose monohydrate	292996	1
Bombastus	Lactose monohydrate	282696/0	1
Bombastus	Lactose monohydrate	284321/0	1
Bombastus	Lactose monohydrate	293074	1
Bombastus	Lactose monohydrate	295714	1
Bombastus	Lactose monohydrate	311280/0	1
Bombastus	Lactose monohydrate	309413/0	1
Bombastus	Lactose monohydrate	286697/0	1
Bombastus	Lactose monohydrate	288015/0	1
Bombastus	Lactose monohydrate	294154	1
Bombastus	Lactose monohydrate	293290	1
Bombastus	Lactose monohydrate	292982	2
Bombastus	Lactose monohydrate	305336	2
Bombastus	Lactose monohydrate	294897/0	1
Bombastus	Lactose monohydrate	305227/0	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Lactose monohydrate	294184/0	1
Bombastus	Lactose monohydrate	295011/0	1
Bombastus	Lactose monohydrate	307573/0	1
Bombastus	Lactose monohydrate	306989	2
Bombastus	Lactose monohydrate	293918/0	1
Bombastus	Lactose monohydrate	310912/0	1
Bombastus	Lactose monohydrate	306964	1
Bombastus	Lactose monohydrate	307444/0	1
Bombastus	Lactose monohydrate	275194/0	1
Bombastus	Lactose monohydrate	307543	2
Bombastus	Lactose monohydrate	307741	1
Bombastus	Lactose monohydrate	272680/0	1
Bombastus	Lactose monohydrate	274926	1
Bombastus	Lactose monohydrate	308117	1
Bombastus	Lactose monohydrate	3077410	1
Bombastus	Lactose monohydrate	306990	1
Bombastus	Lactose monohydrate	306652	1
Bombastus	Lactose monohydrate	306850	1
Bombastus	Lactose monohydrate	275225/0	1
Bombastus	Lactose monohydrate	285087/0	1
Bombastus	Lactose monohydrate	309509	1
Bombastus	Lactose monohydrate	310082/0	1
Bombastus	Lactose monohydrate	308438	1
Bombastus	Lactose monohydrate	17001212	1
Bombastus	Lactose monohydrate	17001199	1
Bombastus	Lactose monohydrate	300657	2
Bombastus	Lactose monohydrate	310905	1
Bombastus	Lactose monohydrate	310691/0	1
Bombastus	Lactose monohydrate	311362/0	1
Bombastus	Lactose monohydrate	308117/0	1
Bombastus	Lactose monohydrate	309351	1
Bombastus	Lactose monohydrate	304755	1
Bombastus	Lactose monohydrate	309508	1
Bombastus	Lactose monohydrate	309509/0	1
Bombastus	Lactose monohydrate	3020165	1
Bombastus	Lactose monohydrate	293289	1
Bombastus	Lactose monohydrate	303658	2
Bombastus	Lactose monohydrate	17000754	1
Bombastus	Lactose monohydrate	293970	1
Bombastus	Lactose monohydrate	311280	1
Bombastus	Lactose monohydrate	310691	1
Bombastus	Lactose monohydrate	310905/0	1
Bombastus	Lactose monohydrate	289238	1
Bombastus	Lactose monohydrate	288927	1
Bombastus	Lactose monohydrate	291464	1
Bombastus	Lactose monohydrate	2932880	1
Bombastus	Lactose monohydrate	293649	1
Bombastus	Lactose monohydrate	295476	1
Bombastus	Lavender flower	268832	1
Bombastus	Lavender flower	279213	1
Bombastus	Lavender flower	274469	1
Bombastus	Lidocaine hydrochlori...	180831-01	1
Bombastus	Lime flower, cut	275129	2
Bombastus	Lime flower, cut	280678	3
Bombastus	Linseed	282734	3
Bombastus	Linseed	274032	1
Bombastus	Linseed	278457	6
Bombastus	Linseed	269500	1
Bombastus	Linseed	1610	1
Bombastus	Linseed	295376	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Liquorice root, peeled	268001	1
Bombastus	Liquorice root, peeled	270192	1
Bombastus	Liquorice root, natur...	903312151	1
Bombastus	Liquorice root, natur...	274921	3
Bombastus	Liquorice root, natur...	5091305	1
Bombastus	Liquorice root, natur...	2901312151	1
Bombastus	Liquorice root, natur...	276468	6
Bombastus	Liquorice root, natur...	3713Q2151	1
Bombastus	Liquorice root, natur...	278377	3
Bombastus	Liquorice root, natur...	285185	2
Bombastus	Liquorice root, natur...	286584	1
Bombastus	Liquorice root, natur...	279409	1
Bombastus	Lovage root, cut	273001	1
Bombastus	Macrogols 4000	13136016	1
Bombastus	Magnesium sulfate hep...	277424	1
Bombastus	Maize style, cut	285912	1
Bombastus	Maize style, cut	282574	1
Bombastus	Maize style, cut	276776	1
Bombastus	Maize style, cut	81113a	1
Bombastus	Male speedwell wort, ...	274238	1
Bombastus	Mallow flower, cut	13267304	1
Bombastus	Mallow flower, cut	273803	1
Bombastus	Marian thistle fruit	278379	1
Bombastus	Marian thistle fruit	279405	3
Bombastus	Marian thistle fruit	277720	3
Bombastus	Marian thistle fruit	296193	1
Bombastus	Marian thistle fruit	280304	1
Bombastus	Marian thistle fruit	282283	1
Bombastus	Marian thistle fruit	13122013a	1
Bombastus	Marian thistle fruit	299483	1
Bombastus	Marian thistle fruit	282034	2
Bombastus	Marian thistle fruit	278378	1
Bombastus	Marian thistle fruit	284420	1
Bombastus	Marian thistle fruit	269499	1
Bombastus	Marian thistle fruit	296535	1
Bombastus	Marshmallow root, nat...	270467	1
Bombastus	Marshmallow root, nat...	276741	1
Bombastus	Marshmallow root, pee...	273458	1
Bombastus	Matricaria flower	287855	1
Bombastus	Matricaria flower	276307	1
Bombastus	Matricaria flower	264997	1
Bombastus	Matricaria flower	281103	2
Bombastus	Matricaria flower	271781	1
Bombastus	Matricaria flower	276095	5
Bombastus	Matricaria flower	277959	2
Bombastus	Matricaria flower	299185	1
Bombastus	Matricaria flower	283710	2
Bombastus	Matricaria flower	268712	1
Bombastus	Matricaria flower	268408	1
Bombastus	Matricaria flower	274106	1
Bombastus	Matricaria flower	265419	1
Bombastus	Matricaria flower	268408	1
Bombastus	Matricaria flower	274105	1
Bombastus	Matricaria flower	273667	1
Bombastus	Melissa leaf, cut	270044	1
Bombastus	Melissa leaf, cut	271499	2
Bombastus	Melissa leaf, cut	277722	1
Bombastus	Methyl-4-hydroxybenzo...	1212011-01	2
Bombastus	Metronidazole	123Validierung	1
Bombastus	Metronidazole	15253032	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Metronidazole	18026401	6
Bombastus	Miconazole nitrate	15I08-B07-318687	3
Bombastus	Miconazole nitrate	16144725	14
Bombastus	Mistletoe, cut	277432	2
Bombastus	Nettle leaf, cut	270987	3
Bombastus	Nettle leaf, cut	273229	2
Bombastus	Nettle leaf, cut	270854	3
Bombastus	Nettle leaf, cut	275848	3
Bombastus	Nettle leaf, cut	275763	1
Bombastus	Nettle leaf, cut	266638	1
Bombastus	Nettle leaf, cut	280305	2
Bombastus	Nettle leaf, cut	276407	2
Bombastus	Nettle leaf, cut	266099	1
Bombastus	Non-ionic emulsifying...	295568	6
Bombastus	Non-ionic emulsifying...	306665	1
Bombastus	Non-ionic emulsifying...	278281	4
Bombastus	Non-ionic emulsifying...	274059	2
Bombastus	Non-ionic emulsifying...	267575	2
Bombastus	Non-ionic emulsifying...	306666	8
Bombastus	Non-ionic emulsifying...	284009	2
Bombastus	Non-ionic emulsifying...	2995568	1
Bombastus	Non-ionic emulsifying...	292319	8
Bombastus	Non-ionic emulsifying...	3066	1
Bombastus	Non-ionic emulsifying...	304997	5
Bombastus	Norethisterone acetate	13J08-N04	2
Bombastus	Oak bark, cut	274392	1
Bombastus	Oak bark, cut	279525	2
Bombastus	Oak bark, cut	272704	1
Bombastus	Oak bark, cut	274611	1
Bombastus	Oak bark, cut	265717	1
Bombastus	Oak bark, cut	229090	1
Bombastus	Oak bark, cut	286297	1
Bombastus	Oxytetracycline hydro...	13L19-B33-290732	1
Bombastus	Paraffin, hard	293848	1
Bombastus	Peppermint leaf, cut	255719	1
Bombastus	Peppermint leaf, cut	279214	1
Bombastus	Peppermint leaf, cut	281282	1
Bombastus	Peppermint leaf, cut	279634	1
Bombastus	Peppermint leaf, cut	273843	2
Bombastus	Peppermint leaf, cut	266772	2
Bombastus	Peppermint leaf, cut	285914	1
Bombastus	Peppermint leaf, cut	279914	2
Bombastus	Peppermint leaf, cut	281043	3
Bombastus	Peppermint leaf, cut	275295	3
Bombastus	Peppermint leaf, cut	246514	1
Bombastus	Peppermint leaf, cut	276514	1
Bombastus	Peppermint leaf, cut	282509	1
Bombastus	Peppermint leaf, cut	272956	3
Bombastus	Peppermint leaf, cut	271499	1
Bombastus	Peppermint leaf, cut	278679	2
Bombastus	Potato starch	12214303	2
Bombastus	Prednisolone acetate	2255	1
Bombastus	Prednisolone, microni...	17F07-B05-341574	5
Bombastus	Prednisolone, microni...	17F07-B05-341578	4
Bombastus	Raspberry leaf, cut	248586	1
Bombastus	Raspberry leaf, cut	272898	1
Bombastus	Raspberry leaf, cut	269870	1
Bombastus	Raspberry leaf, cut	273668	2
Bombastus	Raspberry leaf, cut	279686	3
Bombastus	Raspberry leaf, cut	267122	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Raspberry leaf, cut	273778	1
Bombastus	Raspberry leaf, cut	274507	1
Bombastus	Restharrow root, cut	272831	2
Bombastus	Restharrow root, cut	274421	1
Bombastus	Restharrow root, cut	283708	1
Bombastus	Rhubarb, powdered	11043806	1
Bombastus	Ribwort plantain, cut	276453	1
Bombastus	Ribwort plantain, cut	275479	1
Bombastus	Rose hips peel, cut	283458	1
Bombastus	Rose hips peel, cut	285184	1
Bombastus	Rose hips peel, cut	275394	2
Bombastus	Rose hips peel, cut	231013A	1
Bombastus	Rose hips peel, cut	277431	3
Bombastus	Rose hips peel, cut	257166	1
Bombastus	Rosemary leaf, cut	270171	1
Bombastus	Rosemary leaf, cut	277128	1
Bombastus	Rosemary leaf, cut	281843	1
Bombastus	Rosemary leaf, cut	281248	1
Bombastus	Sage leaf, cut	270115	1
Bombastus	Sage leaf, cut	277512	1
Bombastus	Sage leaf, cut	255381	1
Bombastus	Sage leaf, cut	271253	1
Bombastus	Salicylic acid	12323403	4
Bombastus	Salicylic acid	16050017	12
Bombastus	Senna leaf, cut	273547	3
Bombastus	Senna leaf, cut	273474	1
Bombastus	Senna leaf, cut	278109	1
Bombastus	Senna leaf, cut	280581	2
Bombastus	Senna leaf, cut	280448	1
Bombastus	Shepherd's purse, cut	12345	2
Bombastus	Shepherd's purse, cut	269977	1
Bombastus	Shepherd's purse, cut	284136	1
Bombastus	Shepherd's purse, cut	278856	2
Bombastus	Shepherd's purse, cut	269594	4
Bombastus	Sodium benzoate	17353011	3
Bombastus	Sodium citrate	73155098	1
Bombastus	Sodium cromoglicate	15414502	1
Bombastus	Sodium edetate		1
Bombastus	Sodium monohydrogenph...	11H04-N06	1
Bombastus	Sorbitol	13278206	3
Bombastus	St. John's wort, cut	275106	6
Bombastus	St. John's wort, cut	265761	1
Bombastus	St. John's wort, cut	268834	2
Bombastus	St. John's wort, cut	281530	1
Bombastus	St. John's wort, cut	268913	1
Bombastus	Sucrose	15156008	1
Bombastus	Tall golden rod, cut	271459	1
Bombastus	Tall golden rod, cut	278697	2
Bombastus	Tall golden rod, cut	209974	1
Bombastus	Tall golden rod, cut	274290	1
Bombastus	Tall golden rod, cut	271460	1
Bombastus	Tall golden rod, cut	270047	2
Bombastus	Tall golden rod, cut	261756	1
Bombastus	Thyme, rubbed	243485	1
Bombastus	Thyme, rubbed	273909	3
Bombastus	Thyme, rubbed	277996	2
Bombastus	Thyme, rubbed	277466	1
Bombastus	Tretinoin	1402007-03	4
Bombastus	Triclosan	13L18-35-289984	1
Bombastus	Turmeric, javanese, cut	277130	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Bombastus	Turmeric, javanese, cut	282221	2
Bombastus	Turmeric, javanese, cut	288629	1
Bombastus	Turmeric, javanese, cut	28221	1
Bombastus	Turmeric, javanese, cut	262792	1
Bombastus	Turmeric, javanese, cut	284242	2
Bombastus	Urea	16107424	46
Bombastus	Urea	11150721	2
Bombastus	Valerian root, cut	274722	1
Bombastus	Valerian root, cut	267863	1
Bombastus	Valerian root, cut	269457	1
Bombastus	Walnut bark, cut	711B1206130	1
Bombastus	Wax, yellow	caelo	1
Bombastus	Wax, yellow	279999	1
Bombastus	Wax, yellow	267702	1
Bombastus	Wax, yellow	270682	3
Bombastus	Wax, yellow	285109	1
Bombastus	Wax, yellow	292025	1
Bombastus	Wax, yellow	296261	1
Bombastus	Wild pansy, cut	278106	2
Bombastus	Wild pansy, cut	280890	1
Bombastus	Wild pansy, cut	220232	1
Bombastus	Wild pansy, cut	272832	1
Bombastus	Wild pansy, cut	282661	2
Bombastus	Wild pansy, cut	278106	2
Bombastus	Wild pansy, cut	279314	3
Bombastus	Wild pansy, cut	276259	3
Bombastus	Wild pansy, cut	274237	2
Bombastus	Willow bark, cut	277847	1
Bombastus	Willow bark, cut	261615	1
Bombastus	Willow bark, cut	161594	1
Bombastus	Willow bark, cut	282270	1
Bombastus	Willow bark, cut	288869	1
Bombastus	Willow bark, cut	278458	2
Bombastus	Willow bark, cut	281533	2
Bombastus	Willow bark, cut	277005	3
Bombastus	Willow bark, cut	251385	1
Bombastus	Willowherb, small flo. . .	281846	2
Bombastus	Wormwood herb, cut	281000	1
Bombastus	Wormwood herb, cut	221868	1
Bombastus	Wormwood herb, cut	265417	2
Bombastus	Yarrow herb, cut	15006702	1
Bombastus	Yarrow herb, cut	274284	1
Bombastus	Yarrow herb, cut	276498	1
Bombastus	Yarrow herb, cut	275964	1
Bombastus	Yarrow herb, cut	274610	2
Borax Francais	Sodium tetraborate	1402002-01	1
buchler GmbH	Quinidine sulfate dihydrochloride	16a06-b04-340797	1
Burton Feingold	Calcium ascorbate	L008415	1
C	Robert herb, cut	11147810	3
Caelo	Sodium sulfate decahydrate	11294603	2
Cae	Bean shells, cut	Int.Ch.313011	1
Cae	Bilberry fruit, dried	212172	1
Cae	Bitter fennel	Int.Ch.314090	1
Caelo/Phönix	Glucose monohydrate	15095601	13
Caelo	5-aminolevulinic acid. . .	11355616	1
Caelo	5-aminolevulinic acid. . .	12093205	2
Caelo	5-aminolevulinic acid. . .	12093208	1
Caelo	5-aminolevulinic acid. . .	13017402	1
Caelo	5-aminolevulinic acid. . .	13017403	1
Caelo	5-aminolevulinic acid. . .	14145003	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	5-aminolevulinic acid...	14145007	2
Caelo	5-aminolevulinic acid...	14145011	1
Caelo	5-aminolevulinic acid...	15254803	2
Caelo	Abrotani herb, cut	2091305	1
Caelo	Acetylcysteine	12299004	1
Caelo	Acetylcysteine	10230713	1
Caelo	Acetylcysteine	12299021	1
Caelo	Acetylcysteine	12299027	1
Caelo	Acetylcysteine	130702601	1
Caelo	Acetylcysteine	15114007	2
Caelo	Acetylcysteine	15114013	2
Caelo	Acetylcysteine	16139802	2
Caelo	Acetylcysteine	16139808	1
Caelo	Acetylcysteine	1011171701	1
Caelo	Acetylcysteine	16334006	3
Caelo	Acetylsalicylic acid	15334910	1
Caelo	Acetylsalicylic acid	15334907	4
Caelo	Acetylsalicylic acid	15081308	2
Caelo	Acetylsalicylic acid	15081304	1
Caelo	Acetylsalicylic acid	15122701	1
Caelo	Acetylsalicylic acid	15081310	1
Caelo	Acetylsalicylic acid	301115	1
Caelo	Acetylsalicylic acid	15081301	3
Caelo	Acetylsalicylic acid	15334902	3
Caelo	Acetylsalicylic acid	18122612	1
Caelo	Acetylsalicylic acid	16131607	2
Caelo	Acetylsalicylic acid	1110417	1
Caelo	Acetylsalicylic acid	17023107	7
Caelo	Acetylsalicylic acid	17023108	2
Caelo	Acetylsalicylic acid	1061117/17023107	1
Caelo	Acetylsalicylic acid	17023111	1
Caelo	Acetylsalicylic acid	17023110	6
Caelo	Acetylsalicylic acid	17329204	4
Caelo	Acetylsalicylic acid	18122602	5
Caelo	Acetylsalicylic acid	18108306	1
Caelo	Acetylsalicylic acid	18122605	2
Caelo	Acetylsalicylic acid	17023102	7
Caelo	Acetylsalicylic acid	13324116	1
Caelo	Acetylsalicylic acid	15334917	4
Caelo	Acetylsalicylic acid	12212223	3
Caelo	Acetylsalicylic acid	10099907	1
Caelo	Acetylsalicylic acid	12212219	1
Caelo	Acetylsalicylic acid	11257126	1
Caelo	Acetylsalicylic acid	12212203	1
Caelo	Acetylsalicylic acid	14189603	1
Caelo	Acetylsalicylic acid	12212230	1
Caelo	Acetylsalicylic acid	12212234	2
Caelo	Acetylsalicylic acid	12212233	1
Caelo	Acetylsalicylic acid	12212236	2
Caelo	Acetylsalicylic acid	13324101	1
Caelo	Acetylsalicylic acid	13324109	4
Caelo	Acetylsalicylic acid	13324119	1
Caelo	Acetylsalicylic acid	1703010-02	1
Caelo	Acetylsalicylic acid	13297511	1
Caelo	Acetylsalicylic acid	13324118	1
Caelo	Acetylsalicylic acid	14250101	2
Caelo	Acetylsalicylic acid	13324105	1
Caelo	Aciclovir	18118-B04-355420	1
Caelo	Acorus root, natural...	12066501	1
Caelo	Acorus root, natural...	13127801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Acorus root, natural,...	284875	2
Caelo	Acorus root, peeled	313060	1
Caelo	Acorus root, peeled	12295903	1
Caelo	Acorus root, peeled	28121205	1
Caelo	Acorus root, peeled	24081203	1
Caelo	Acorus root, peeled	15061202	1
Caelo	Acorus root, peeled	12056701	1
Caelo	Acorus root, peeled	241013010	1
Caelo	Acorus root, peeled	12295902	4
Caelo	Acorus root, peeled	12056703	1
Caelo	Acorus root, peeled	8031306	1
Caelo	Acorus root, peeled	12066501	7
Caelo	Acorus root, peeled	12066502	1
Caelo	Acorus root, peeled	12275604-3544411013	1
Caelo	Acorus root, peeled	13101306	1
Caelo	Acorus root, peeled	12275604	2
Caelo	Acorus root, peeled	12239008	1
Caelo	Acorus root, peeled	12295901	1
Caelo	Acorus root, peeled	12295904	4
Caelo	Acorus root, peeled	7061302	1
Caelo	Acorus root, peeled	12275602	7
Caelo	Acorus root, peeled	6701B130103	21
Caelo	Agarikon	11024502	1
Caelo	Agnus castus fruit	13193404	1
Caelo	Agnus castus fruit	12274009	3
Caelo	Agnus castus fruit	12016905	2
Caelo	Agnus castus fruit	3413A-03015	1
Caelo	Agnus castus fruit	15065105	1
Caelo	Agnus castus fruit	15318705	1
Caelo	Agnus castus fruit	80175158	2
Caelo	Agnus castus fruit	13193406	1
Caelo	Agnus castus fruit	15318718	1
Caelo	Agnus castus fruit	12274014	1
Caelo	Agnus castus fruit	12274013	2
Caelo	Agnus castus fruit	13193401	2
Caelo	Agnus castus fruit	13193410	2
Caelo	Agnus castus fruit	13193403	1
Caelo	Agnus castus fruit	14194408	1
Caelo	Agnus castus fruit	15318702	1
Caelo	Agnus castus fruit	1513546	1
Caelo	Agnus castus fruit	12274012	1
Caelo	Agnus castus fruit	12274008	1
Caelo	Agnus castus fruit	12274010	10
Caelo	Agnus castus fruit	12016908	1
Caelo	Agnus castus fruit	12274002	1
Caelo	Agnus castus fruit	12274003	4
Caelo	Agnus castus fruit	11026318	1
Caelo	Agnus castus fruit	12016902	1
Caelo	Agnus castus fruit	12016909	4
Caelo	Agnus castus fruit	11026309	1
Caelo	Agrimony, cut	11374005	2
Caelo	Agrimony, cut	12382901	3
Caelo	Agrimony, cut	12382903	1
Caelo	Agrimony, cut	1137407	1
Caelo	Agrimony, cut	10325502	2
Caelo	Agrimony, cut	3481-A-110324	1
Caelo	Agrimony, cut	11374006	6
Caelo	Agrimony, cut	10325505	1
Caelo	Alchemilla alpina her...	12147206	3
Caelo	Alchemilla alpina her...	12147202	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Alchemilla alpina her. ...	11063608	1
Caelo	Alchemilla alpina her. ...	12147203	1
Caelo	Alchemilla alpina her. ...	12147203-3470312131	1
Caelo	Alchemilla alpina her. ...	1214204	1
Caelo	Alchemilla alpina her. ...	12147205	2
Caelo	Alchemilla, cut	11182804	3
Caelo	Alchemilla, cut	3613M01616	1
Caelo	Alchemilla, cut	13141103-3400310161	1
Caelo	Alchemilla, cut	1406211616	1
Caelo	Alchemilla, cut	13141103-3382319003	1
Caelo	Alchemilla, cut	12153306	5
Caelo	Alchemilla, cut	13091304	1
Caelo	Alchemilla, cut	12342901	1
Caelo	Alchemilla, cut	12153307	3
Caelo	Alchemilla, cut	2081304	1
Caelo	Alchemilla, cut	13141102	6
Caelo	Alchemilla, cut	1215309	1
Caelo	Alchemilla, cut	13141103	36
Caelo	Alchemilla, cut	125306	1
Caelo	Alchemilla, cut	13141101	7
Caelo	Alchemilla, cut	5041301	1
Caelo	Alchemilla, cut	H-040513-2	1
Caelo	Alchemilla, cut	11182808	1
Caelo	Alchemilla, cut	intch:318080	1
Caelo	Alchemilla, cut	11063605	1
Caelo	Alchemilla, cut	13141105	3
Caelo	Alchemilla, cut	1213309	1
Caelo	Alchemilla, cut	11182805	1
Caelo	Alchemilla, cut	1118206	1
Caelo	Alchemilla, cut	21081203	1
Caelo	Alchemilla, cut	1100912	1
Caelo	Alchemilla, cut	11002906	5
Caelo	Alchemilla, cut	1215305	2
Caelo	Alchemilla, cut	12153301	3
Caelo	Alchemilla, cut	15147002	2
Caelo	Alchemilla, cut	11182807	2
Caelo	Alchemilla, cut	8121201	1
Caelo	Alchemilla, cut	14066008	1
Caelo	Alchemilla, cut	11182801	3
Caelo	Alchemilla, cut	14066003	1
Caelo	Alchemilla, cut	131411002	1
Caelo	Alchemilla, cut	13141103-3456311103	1
Caelo	Alchemilla, cut	12153303	7
Caelo	Alchemilla, cut	5013I1616	1
Caelo	Alchemilla, cut	9011405	1
Caelo	Alchemilla, cut	13141106	8
Caelo	Alchemilla, cut	27011401	1
Caelo	Alchemilla, cut	12153305	25
Caelo	Alchemilla, cut	13141106-3545412081	1
Caelo	Alchemilla, cut	14066002	2
Caelo	Alfatradiol	16087511	1
Caelo	Alfatradiol	16180409	3
Caelo	Alfatradiol	160875105728	1
Caelo	Alfatradiol	14318511	1
Caelo	Alfatradiol	16087512	4
Caelo	Alfatradiol	17G03-B11-341398	3
Caelo	Alfatradiol	16087504	8
Caelo	Alfatradiol	16087507	1
Caelo	Alfatradiol	4318512	1
Caelo	Alfatradiol	19A24-B03-360709	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Alfatradiol	14318512	3
Caelo	Alfatradiol	16180408	3
Caelo	Alfatradiol	14318509	1
Caelo	Alfatradiol	14139103	5
Caelo	Alfatradiol	14037904	3
Caelo	Alfatradiol	13130808	1
Caelo	Alfatradiol	511211	1
Caelo	Alfatradiol	13130806	2
Caelo	Alfatradiol	13130805	5
Caelo	Alfatradiol	111161607	1
Caelo	Alfatradiol	12164207	7
Caelo	Alfatradiol	14139104	11
Caelo	Alfatradiol	10319904	1
Caelo	Alfatradiol	12H29-N05	6
Caelo	Alfatradiol	10319909	1
Caelo	Alfatradiol	12164201	3
Caelo	Alfatradiol	15356804	7
Caelo	Alfatradiol	11061004	1
Caelo	Alfatradiol	11061005	1
Caelo	Alfatradiol	11161805	4
Caelo	Alfatradiol	11061001	3
Caelo	Alfatradiol	11161807	3
Caelo	Alfatradiol	14C20-B02-294944	2
Caelo	Alfatradiol	5051406	1
Caelo	Alfatradiol	14139104	2
Caelo	Alfatradiol	14318502	1
Caelo	Alfatradiol	14139102	1
Caelo	Alfatradiol	15356801	16
Caelo	Alfatradiol	14139102	1
Caelo	Alfatradiol	5107516	1
Caelo	Alfatradiol	14139103	1
Caelo	Alfatradiol	14139105	2
Caelo	Alfatradiol	15380001	1
Caelo	Alfatradiol	14I22-B30-302019	1
Caelo	Alfatradiol	15380001	7
Caelo	Alfatradiol	14130104	1
Caelo	Alfatradiol	2061504	1
Caelo	Alfatradiol	1536801	1
Caelo	Allantoin	16228906	11
Caelo	Allantoin	16228904	1
Caelo	Allantoin	16228915	21
Caelo	Allantoin	16228914	1
Caelo	Allantoin	17355903	4
Caelo	Allantoin	18158602	1
Caelo	Allantoin	81542358	2
Caelo	Allantoin	15220006	15
Caelo	Allantoin	14247609	7
Caelo	Allantoin	121127	1
Caelo	Allantoin	15220010	3
Caelo	Allantoin	15220012	10
Caelo	Allantoin	15220011	2
Caelo	Allantoin	15275302	1
Caelo	Allantoin	152200012	1
Caelo	Allantoin	15275305	2
Caelo	Allantoin	15220007	4
Caelo	Allantoin	15275312	8
Caelo	Allantoin	16228903	9
Caelo	Allantoin	15275304	2
Caelo	Allantoin	12374418	2
Caelo	Allantoin	1522000	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Allantoin	15220004	1
Caelo	Allantoin	28101702	1
Caelo	Allantoin	11247504	5
Caelo	Allantoin	20071202	1
Caelo	Allantoin	15220005	2
Caelo	Allantoin	1522008	1
Caelo	Allantoin	12374423	2
Caelo	Allantoin	12251510	12
Caelo	Allantoin	140512	1
Caelo	Allantoin	12251511	1
Caelo	Allantoin	10318305	1
Caelo	Allantoin	10064012	1
Caelo	Allantoin	10064008	4
Caelo	Allantoin	2264594	1
Caelo	Allantoin	12374401	1
Caelo	Allantoin	12374404	2
Caelo	Allantoin	10064015	1
Caelo	Allantoin	12374407	4
Caelo	Allantoin	2713448	1
Caelo	Allantoin	4247602	1
Caelo	Allantoin	5051402	1
Caelo	Allantoin	1251505	1
Caelo	Allantoin	11247510	2
Caelo	Allantoin	12374414	8
Caelo	Allantoin	12374408	1
Caelo	Allantoin	10318312	1
Caelo	Allantoin	12374422	9
Caelo	Allantoin	10318313	1
Caelo	Allantoin	1113E-01087	1
Caelo	Allantoin	2713431	1
Caelo	Allantoin	1225101	1
Caelo	Allantoin	12374417	1
Caelo	Allantoin	12251505	9
Caelo	Allantoin	12251502	2
Caelo	Allantoin	12251501	5
Caelo	Allantoin	11247511	2
Caelo	Allantoin	20140612-15	1
Caelo	Allantoin	12374415	12
Caelo	Allantoin	15220008	10
Caelo	Almonds, bitter, whole	13242001	1
Caelo	Almonds, bitter, whole	16221202	1
Caelo	Almonds, bitter, whole	15193103	1
Caelo	Almonds, bitter, whole	5121301	1
Caelo	Almonds, bitter, whole	12239604	1
Caelo	Aloes, cape, cut	15283001	1
Caelo	Aloes, cape, cut	11061406	1
Caelo	Aloes, cape, cut	15283010	2
Caelo	Aloes, cape, cut	16202401	1
Caelo	Aloes, cape, cut	16202404	1
Caelo	Aloes, cape, cut	11041402	1
Caelo	Aloes, cape, cut	16202403	1
Caelo	Aloes, cape, cut	25111303	1
Caelo	Aloes, cape, cut	12338206	2
Caelo	Aloes, cape, cut	505125	1
Caelo	Aloes, cape, cut	12232102	1
Caelo	Aloes, cape, cut	12338201	4
Caelo	Aloes, cape, cut	12338203	4
Caelo	Aloes, cape, cut	12338209	1
Caelo	Aloes, cape, cut	25071402	1
Caelo	Aloes, cape, powdered	16202410	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Aloes, cape, powdered	15206514	1
Caelo	Aloes, cape, powdered	15206516	1
Caelo	Aloes, cape, powdered	15206507	1
Caelo	Aloes, cape, powdered	14065905	1
Caelo	Aloes, cape, powdered	14065906	1
Caelo	Aloes, cape, powdered	12183013	3
Caelo	Aloes, cape, powdered	12183008	1
Caelo	Aloes, cape, powdered	12183009	3
Caelo	Aloes, cape, powdered	11132116	2
Caelo	Aloes, cape, powdered	11132116	1
Caelo	Aloes, cape, powdered	12183001	1
Caelo	Aloes, cape, powdered	11132115	1
Caelo	Aloes, cape, powdered	12183007	1
Caelo	Aloes, cape, powdered	15206501	1
Caelo	alpha-tocopheryl acet...	1705003-01	1
Caelo	Aluminium chloride he...	12114912	4
Caelo	Aluminium chloride he...	3113E-01095	1
Caelo	Aluminium chloride he...	13G03-N02	21
Caelo	Aluminium chloride he...	5081301	1
Caelo	Aluminium chloride he...	28121301	1
Caelo	Aluminium chloride he...	12114918	17
Caelo	Aluminium chloride he...	12323710	8
Caelo	Aluminium chloride he...	12323711	17
Caelo	Aluminium chloride he...	8031301	1
Caelo	Aluminium chloride he...		2
Caelo	Aluminium chloride he...	19041302	1
Caelo	Aluminium chloride he...	10251314	1
Caelo	Aluminium chloride he...	9071303	1
Caelo	Aluminium chloride he...	18061308	1
Caelo	Aluminium chloride he...	12127-N01	1
Caelo	Aluminium chloride he...	51/AD1	1
Caelo	Aluminium chloride he...	15051302	1
Caelo	Aluminium chloride he...	12114920	11
Caelo	Aluminium chloride he...	12021303	1
Caelo	Aluminium chloride he...	11229907	1
Caelo	Aluminium chloride he...	12114910	6
Caelo	Aluminium chloride he...	11229917	3
Caelo	Aluminium chloride he...	11229914	1
Caelo	Aluminium chloride he...	11229905	1
Caelo	Aluminium chloride he...	12047503	2
Caelo	Aluminium chloride he...	12114922	14
Caelo	Aluminium chloride he...	15080715	1
Caelo	Aluminium chloride he...	13L13-B03-287298	1
Caelo	Aluminium chloride he...	12047514	3
Caelo	Aluminium chloride he...	10251307	1
Caelo	Aluminium chloride he...	15080701	1
Caelo	Aluminium chloride he...	3091205	1
Caelo	Aluminium chloride he...	12114908	8
Caelo	Aluminium chloride he...	12114906	6
Caelo	Aluminium chloride he...	12323717	5
Caelo	Aluminium chloride he...	12114904	2
Caelo	Aluminium chloride he...	23111203	1
Caelo	Aluminium chloride he...	11229906	1
Caelo	Aluminium chloride he...	13/13-B03-287298	1
Caelo	Aluminium chloride he...	12114919	1
Caelo	Aluminium potassium s...	12022401	1
Caelo	Aluminium potassium s...	12184408	3
Caelo	Aluminium potassium s...	12184406	2
Caelo	Aluminium potassium s...	11218611	2
Caelo	Aluminium potassium s...	12144607	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Aluminium potassium s...	12022406	1
Caelo	Aluminium potassium s...	12239102	2
Caelo	Aluminium potassium s...	12184401	1
Caelo	Aluminium potassium s...	1218411	1
Caelo	Aluminium potassium s...	13025401	3
Caelo	Aluminium potassium s...	12184411	1
Caelo	Aluminium potassium s...	11254805	1
Caelo	Aluminium potassium s...	12144605	4
Caelo	Aluminium potassium s...	13025404	3
Caelo	Aluminium potassium s...	13098702	1
Caelo	Aluminium potassium s...	14263804	1
Caelo	Aluminium potassium s...	15172501	1
Caelo	Aluminium potassium s...	90033229	1
Caelo	Aluminium potassium s...	2713Q1092	1
Caelo	Aluminium potassium s...	12 24-M03	1
Caelo	Aluminium potassium s...	12022444	1
Caelo	Aluminium potassium s...	12184401	1
Caelo	Aluminium potassium s...	12144602	1
Caelo	Aluminium potassium s...	12184407	3
Caelo	Aluminium potassium s...	12022409	1
Caelo	Aluminium potassium s...	12184402	3
Caelo	Aluminium potassium s...	12184407	5
Caelo	Aluminium potassium s...	12022406	2
Caelo	Aluminium potassium s...	12144603	1
Caelo	Aluminium potassium s...	12022404	1
Caelo	Aluminium potassium s...	11017211	1
Caelo	Aluminium potassium s...	510211092	1
Caelo	Aluminium potassium s...	504311092	1
Caelo	Ambroxol hydrochloride	18013208	1
Caelo	Ambroxol hydrochloride	18013211	4
Caelo	Ambroxol hydrochloride	1801306	1
Caelo	Ambroxol hydrochloride	17158204	1
Caelo	Ambroxol hydrochloride	18013205	5
Caelo	Ambroxol hydrochloride	17002904	1
Caelo	Ambroxol hydrochloride	19061804	3
Caelo	Ambroxol hydrochloride	17286001	2
Caelo	Ambroxol hydrochloride	27101707	1
Caelo	Ambroxol hydrochloride	17158202	2
Caelo	Ambroxol hydrochloride	170522AH/16309201	1
Caelo	Ambroxol hydrochloride	16309201	1
Caelo	Ambroxol hydrochloride	16199602	1
Caelo	Ambroxol hydrochloride	16199603	2
Caelo	Ambroxol hydrochloride	17158201	2
Caelo	Ambroxol hydrochloride	19061803	1
Caelo	Ambroxol hydrochloride	17002902	1
Caelo	Amifampridine	14H11-B06-298447	2
Caelo	Amifampridine	13f05-n01	1
Caelo	Ammonium chloride	12356101	3
Caelo	Ammonium chloride	12240302	4
Caelo	Ammonium chloride	10176319	2
Caelo	Ammonium chloride	10176309	1
Caelo	Ammonium hydrogen car...	25111308	1
Caelo	Ammonium hydrogen car...	11145612	2
Caelo	Ammonium hydrogen car...	9139402	1
Caelo	Ammonium hydrogen car...	12022308	2
Caelo	Ammonium hydrogen car...	12322201	1
Caelo	Ammonium hydrogen car...	13104601	1
Caelo	Ammonium hydrogen car...	12322204	1
Caelo	Ammonium hydrogen car...	23011404	1
Caelo	Ammonium hydrogen car...	72715338	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ammonium hydrogen car...	13104611	1
Caelo	Ammonium hydrogen car...	24101202	1
Caelo	Ammonium hydrogen car...	3712Q-01103	1
Caelo	Ammonium hydrogen car...	3091204	1
Caelo	Angelica root	Int.Ch:318090	1
Caelo	Angelica root	14179506	2
Caelo	Angelica root	10071306	1
Caelo	Angelica root	6081B130109	7
Caelo	Angelica root	269740	1
Caelo	Angelica root	30011302	1
Caelo	Angelica root	12168705	2
Caelo	Angelica root	12350703	9
Caelo	Angelica root	9030301	1
Caelo	Angelica root	5091303	1
Caelo	Angelica root	15101304	1
Caelo	Angelica root	14178506	1
Caelo	Angelica root	12350704	2
Caelo	Angelica root	21011402	1
Caelo	Angelica root	28101402	1
Caelo	Angelica root	12350705	1
Caelo	Angelica root	12080412	1
Caelo	Angelica root	285732	2
Caelo	Angelica root	15162604	2
Caelo	Angelica root	13041711	1
Caelo	Angelica root	12080403	1
Caelo	Angelica root	303299	1
Caelo	Angelica root	16128006	1
Caelo	Angelica root	17053304	1
Caelo	Angelica root	16128002	2
Caelo	Angelica root	12168703	2
Caelo	Angelica root	5162606	1
Caelo	Angelica root	12168701	3
Caelo	Angelica root	17053301	1
Caelo	Angelica root	17053302	1
Caelo	Angelica root	11093304	1
Caelo	Aniseed	11370308	3
Caelo	Aniseed	11121705	1
Caelo	Aniseed	17045601	1
Caelo	Aniseed	13037808-3560412014	1
Caelo	Aniseed	5091304	1
Caelo	Aniseed	2916E-01530	1
Caelo	Aniseed	16165302	1
Caelo	Aniseed	13037805	11
Caelo	Aniseed	11370302	4
Caelo	Aniseed	13037801	1
Caelo	Aniseed	13037804	7
Caelo	Aniseed	13037901	1
Caelo	Aniseed	140913a	1
Caelo	Aniseed	15398607	1
Caelo	Aniseed	13037803	2
Caelo	Aniseed	17045603	1
Caelo	Aniseed	16165306	3
Caelo	Aniseed	16165308	1
Caelo	Aniseed	4061301	1
Caelo	Aniseed	11101202	3
Caelo	Aniseed	111101213	1
Caelo	Aniseed	13037802	9
Caelo	Aniseed	11101203	6
Caelo	Aniseed	13333901	6
Caelo	Aniseed	13037807	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Aniseed	13037809	12
Caelo	Aniseed	5061301	1
Caelo	Aniseed	17045608	1
Caelo	Aniseed	17045611	2
Caelo	Aniseed	12360801	1
Caelo	Aniseed	17045609	1
Caelo	Aniseed	15046103	3
Caelo	Aniseed	13333903	3
Caelo	Aniseed	3314E-01530	1
Caelo	Aniseed	13333907	4
Caelo	Aniseed	5121401	1
Caelo	Aniseed	13333906	3
Caelo	Aniseed	13333909	2
Caelo	Aniseed	13333912	3
Caelo	Aniseed	PrüfNr4913E-01530	1
Caelo	Aniseed	13333910	4
Caelo	Aniseed	1137033	1
Caelo	Aniseed	11101201	1
Caelo	Aniseed	10081207	1
Caelo	Aniseed	11101213	2
Caelo	Aniseed	13333914	1
Caelo	Aniseed	15398602	5
Caelo	Aniseed	15046102	8
Caelo	Aniseed	13333913	1
Caelo	Aniseed	15283201	3
Caelo	Aniseed	15046106	1
Caelo	Aniseed	15046109	5
Caelo	Aniseed	15046101	5
Caelo	Aniseed	15398603	4
Caelo	Aniseed	2215M-01530	1
Caelo	Aniseed	13037808	9
Caelo	Aniseed	11370301	10
Caelo	Aniseed	13333916	1
Caelo	Aniseed	16165304	4
Caelo	Aniseed		1
Caelo	Aniseed	16165305	5
Caelo	Aniseed	11101215	3
Caelo	Aniseed	15046110	1
Caelo	Aniseed	15398605	5
Caelo	Aniseed	13333902	9
Caelo	Aniseed	15049103	1
Caelo	Aniseed	10247202	1
Caelo	Aniseed	16165303	1
Caelo	Aniseed	13333908	7
Caelo	Aniseed	11370303	11
Caelo	Aniseed	11101214	6
Caelo	Aniseed	12226402	5
Caelo	Aniseed, pounded	23041303	1
Caelo	Aniseed, pounded	12111304	1
Caelo	Aniseed, pounded	12360803	2
Caelo	Aniseed, pounded	4041302	1
Caelo	Aniseed, pounded	1314990	1
Caelo	Aniseed, pounded	11242003	1
Caelo	Aniseed, pounded	12111303	3
Caelo	Aniseed, pounded	13149902	1
Caelo	Aniseed, pounded	15149702	1
Caelo	Apple, cut	12227301	1
Caelo	Arginine	14124-B08-303720	1
Caelo	Arginine	15C10-B03-318785	1
Caelo	Arnica flower	12227202	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Arnica flower	12227202	1
Caelo	Arnica flower	11374105	1
Caelo	Arnica flower	11374103	1
Caelo	Arnica flower	1222702	1
Caelo	Arnica flower	IntCh:212032	1
Caelo	Arnica flower	12227203	1
Caelo	Arnica flower	11259804	1
Caelo	Arnica flower	6112013A	1
Caelo	Arnica flower	11285703	1
Caelo	Arnica flower	28071201	1
Caelo	Arnica flower	11374102	2
Caelo	Arnica flower	11259803	1
Caelo	Arnica flower	10207405	1
Caelo	Arnica flower	11285705	1
Caelo	Arnica flower	11374101	2
Caelo	Arnica flower	10081204	1
Caelo	Artichoke leaf, cut	13037607	4
Caelo	Artichoke leaf, cut	13310706	1
Caelo	Artichoke leaf, cut	12173307-3260315072	1
Caelo	Artichoke leaf, cut	13037608	1
Caelo	Artichoke leaf, cut	81232298	1
Caelo	Artichoke leaf, cut	12173306	2
Caelo	Artichoke leaf, cut	13037603	1
Caelo	Artichoke leaf, cut	2312I-03593	1
Caelo	Artichoke leaf, cut	12173307	6
Caelo	Artichoke leaf, cut	12173302	3
Caelo	Artichoke leaf, cut	12173304	3
Caelo	Ascorbic acid	15084503	1
Caelo	Ascorbic acid	13341803/3	1
Caelo	Ascorbic acid	14123915	1
Caelo	Ascorbic acid	14214403	1
Caelo	Ascorbic acid	14123910	1
Caelo	Ascorbic acid	14044701/3	1
Caelo	Ascorbic acid	1071404	1
Caelo	Ascorbic acid	15300002	3
Caelo	Ascorbic acid	12350011	1
Caelo	Ascorbic acid	12272601/1	1
Caelo	Ascorbic acid	22101302	1
Caelo	Ascorbic acid	12350005	2
Caelo	Ascorbic acid	12350003	1
Caelo	Ascorbic acid	12349909	1
Caelo	Ascorbic acid	14123809	2
Caelo	Ascorbic acid	15084507	1
Caelo	Ascorbic acid	16050403	4
Caelo	Ascorbic acid	153000	1
Caelo	Ascorbic acid	18127805	3
Caelo	Ascorbic acid	17360308	2
Caelo	Ascorbic acid	17233201	1
Caelo	Ascorbic acid	17360301	7
Caelo	Ascorbic acid	1234567	1
Caelo	Ascorbic acid	123456	1
Caelo	Ascorbic acid	15359003	2
Caelo	Ascorbic acid	123test	1
Caelo	Ascorbic acid	16345204	1
Caelo	Ascorbic acid	12456	1
Caelo	Ascorbic acid	16345205	6
Caelo	Ascorbic acid	123Test	1
Caelo	Ascorbic acid	16345202	4
Caelo	Ascorbic acid	11208606	2
Caelo	Ascorbic acid	4116E-01016	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ascorbic acid	2061203	1
Caelo	Ascorbic acid	17135803	2
Caelo	Ascorbyl palmitate	16B-B04-525283	1
Caelo	Ascorbyl palmitate	16B10-B04324111	1
Caelo	Ascorbyl palmitate	18B12-B09-525721	3
Caelo	Ascorbyl palmitate	16810-B04525199	1
Caelo	Ascorbyl palmitate	16130-B03-525358	2
Caelo	Ash leaf, cut	12225504	1
Caelo	Ash leaf, cut	12225505	2
Caelo	Ash leaf, cut	12225502	1
Caelo	Asparagus root, cut	13193302	1
Caelo	Asparagus root, cut	13193308	1
Caelo	Atropine sulfate	13357405	1
Caelo	Atropine sulfate	14336807	1
Caelo	Atropine sulfate	16c15b03-319166	1
Caelo	Atropine sulfate	13173809	1
Caelo	Atropine sulfate	15115201	3
Caelo	Atropine sulfate	13357403	1
Caelo	Bacitracin	15C12-B06-317782	2
Caelo	Bacitracin	14A21-B09-292908	1
Caelo	Bacitracin	16I12-B03-331152	1
Caelo	Barberries	12217116	1
Caelo	Barberries	16224202	1
Caelo	Barberries	12217111	1
Caelo	Barberries	12217106	1
Caelo	Barberries	12217109	4
Caelo	Barberries	1221707	1
Caelo	Barberries	12217107	1
Caelo	Barberries	12217103	3
Caelo	Barberries	11058320	1
Caelo	Barberries	12217113	2
Caelo	Barberry root bark, cut	11083509	1
Caelo	Barberry root bark, cut	12112606	5
Caelo	Barberry root bark, cut	12112608	1
Caelo	Barberry root bark, cut	55B110128	1
Caelo	Barberry root bark, cut	12112603	2
Caelo	Barberry root bark, cut	12112605	1
Caelo	Basil, rubbed	13048605	1
Caelo	Basil, rubbed	12007010	4
Caelo	Basil, rubbed	13048601	2
Caelo	Basil, rubbed	13048605-3519411022	1
Caelo	Basil, rubbed	12007008	2
Caelo	Basil, rubbed	12066404	1
Caelo	Bean shells, cut	1103061	1
Caelo	Bean shells, cut	13333803	1
Caelo	Bean shells, cut	1333802	1
Caelo	Bean shells, cut	13333801	1
Caelo	Bean shells, cut	11233603	1
Caelo	Bean shells, cut	29091201	1
Caelo	Bean shells, cut	12358703	5
Caelo	Bean shells, cut	10006507	1
Caelo	Bean shells, cut	11030604	6
Caelo	Bean shells, cut	11030601	1
Caelo	Bean shells, cut	11030605	1
Caelo	Bean shells, cut	6031301	1
Caelo	Bean shells, cut	11030602	5
Caelo	Bean shells, cut	8301004	1
Caelo	Bean shells, cut	14080002	1
Caelo	Bean shells, cut	15060905	1
Caelo	Bearberry leaf, cut	13067805-3390311080	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Bearberry leaf, cut	412A-01518	1
Caelo	Bearberry leaf, cut	20111210	1
Caelo	Bearberry leaf, cut	280889	3
Caelo	Bearberry leaf, cut	4021407	1
Caelo	Bearberry leaf, cut	15051303	1
Caelo	Bearberry leaf, cut	11213011	3
Caelo	Bearberry leaf, cut	11213010	1
Caelo	Bearberry leaf, cut	11213005	4
Caelo	Bearberry leaf, cut	12186502	1
Caelo	Bearberry leaf, cut	15116201	1
Caelo	Bearberry leaf, cut	1218602	1
Caelo	Bearberry leaf, cut	12186504	1
Caelo	Bearberry leaf, cut	5111301	1
Caelo	Bearberry leaf, cut	12186506	3
Caelo	Bearberry leaf, cut	12186503	1
Caelo	Bearberry leaf, cut	4091301	1
Caelo	Bearberry leaf, cut	13067805	6
Caelo	Bearberry leaf, cut	13067802	3
Caelo	Bearberry leaf, cut	17121306	1
Caelo	Bentonite	1307016-01	1
Caelo	Bentonite	1311014-01	1
Caelo	Bentonite	fürHerrSalloch	1
Caelo	Benzocaine	2412M-01118	1
Caelo	Benzocaine	16174008	1
Caelo	Benzocaine	1112706/Nachprüfung	1
Caelo	Benzocaine	15227302	9
Caelo	Benzocaine	16174014	8
Caelo	Benzocaine	16174007	7
Caelo	Benzocaine	16030403	1
Caelo	Benzocaine	11031703	1
Caelo	Benzocaine	16174001	7
Caelo	Benzocaine	13213502	7
Caelo	Benzocaine	15227317	1
Caelo	Benzocaine	3784947	1
Caelo	Benzocaine	13213509	6
Caelo	Benzocaine	15227315	4
Caelo	Benzocaine	15227312	2
Caelo	Benzocaine	15227314	4
Caelo	Benzocaine	15227311	2
Caelo	Benzocaine	15227301	2
Caelo	Benzocaine	15227307	5
Caelo	Benzocaine	1515I-01118	1
Caelo	Benzocaine	13213518	1
Caelo	Benzocaine	13213516	1
Caelo	Benzocaine	14277502	4
Caelo	Benzocaine	13213515	3
Caelo	Benzocaine	5021401	1
Caelo	Benzocaine	6031501	1
Caelo	Benzocaine	161740075429	1
Caelo	Benzocaine	16174020	1
Caelo	Benzocaine	12341906	1
Caelo	Benzocaine	1(19022001)	1
Caelo	Benzocaine	18277801	1
Caelo	Benzocaine	18277805	1
Caelo	Benzocaine	18277806	1
Caelo	Benzocaine	18277802	2
Caelo	Benzocaine	18018804	5
Caelo	Benzocaine	11127506	1
Caelo	Benzocaine	18018806	1
Caelo	Benzocaine	16174019	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Benzocaine	11127508	1
Caelo	Benzocaine	12115007	4
Caelo	Benzocaine	91062169	1
Caelo	Benzocaine	18018808	1
Caelo	Benzocaine	161740195459	1
Caelo	Benzocaine	16174015	1
Caelo	Benzocaine	18041803	1
Caelo	Benzocaine	10072603	3
Caelo	Benzocaine	18018801	7
Caelo	Benzocaine	12341904	7
Caelo	Benzocaine	12115008	2
Caelo	Benzocaine	12115011	3
Caelo	Benzocaine	19031301	1
Caelo	Benzocaine	16174017	1
Caelo	Benzoic acid	16001714	3
Caelo	Benzoic acid	10156104	1
Caelo	Benzoic acid	13074610	4
Caelo	Benzoic acid	13391314	2
Caelo	Benzoic acid	13391320	1
Caelo	Benzoic acid	1339120	1
Caelo	Benzoic acid	15084701	1
Caelo	Benzoic acid	1732701	1
Caelo	Benzoic acid	15084607	4
Caelo	Benzoic acid	15084614	1
Caelo	Benzoic acid	17327101	6
Caelo	Benzoic acid	13391310	1
Caelo	Benzoic acid	15084603	3
Caelo	Betaine monohydrate	4312172831	1
Caelo	Betaine monohydrate	14A14-B02-292983	1
Caelo	Betamethasone dipropi...	17105906	6
Caelo	Betamethasone dipropi...	12114501	1
Caelo	Betamethasone dipropi...	13372503	4
Caelo	Betamethasone dipropi...	13023005	7
Caelo	Betamethasone dipropi...	13372514	2
Caelo	Betamethasone dipropi...	13372515	3
Caelo	Betamethasone dipropi...	13372516	1
Caelo	Betamethasone dipropi...	13236806	1
Caelo	Betamethasone dipropi...	8K03-N13	1
Caelo	Betamethasone dipropi...	13372504	1
Caelo	Betamethasone dipropi...	17353601	11
Caelo	Betamethasone dipropi...	13372516	3
Caelo	Betamethasone dipropi...	15163302	4
Caelo	Betamethasone dipropi...	13372520	1
Caelo	Betamethasone dipropi...	17353603	5
Caelo	Betamethasone dipropi...	15163302	2
Caelo	Betamethasone dipropi...	15163303	9
Caelo	Betamethasone dipropi...	15163301	4
Caelo	Betamethasone dipropi...	2109DM0B00111322	1
Caelo	Betamethasone dipropi...	13372512	1
Caelo	Betamethasone dipropi...	15163301	10
Caelo	Betamethasone dipropi...	15163303	1
Caelo	Betamethasone dipropi...	13372503	3
Caelo	Betamethasone dipropi...	13023002NP	1
Caelo	Betamethasone dipropi...		1
Caelo	Betamethasone dipropi...	1325702	1
Caelo	Betamethasone dipropi...	13125701	4
Caelo	Betamethasone dipropi...	13023001	3
Caelo	Betamethasone dipropi...	13125703	5
Caelo	Betamethasone dipropi...	13236801	6
Caelo	Betamethasone dipropi...	17015905	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Betamethasone dipropi...	1210020-01	1
Caelo	Betamethasone dipropi...	13372504	4
Caelo	Betamethasone dipropi...	12114512	1
Caelo	Betamethasone dipropi...	13236805	4
Caelo	Betamethasone dipropi...	12114507	2
Caelo	Betamethasone dipropi...	17105905	3
Caelo	Betamethasone dipropi...	13236801ok	1
Caelo	Betamethasone dipropi...	12114511	2
Caelo	Betamethasone dipropi...	140507	1
Caelo	Betamethasone dipropi...	12114506	1
Caelo	Betamethasone dipropi...	17105909	2
Caelo	Betamethasone dipropi...	AKD11	1
Caelo	Betamethasone dipropi...	15163204	1
Caelo	Betamethasone dipropi...	13372514	8
Caelo	Betamethasone dipropi...	18090603	7
Caelo	Betamethasone dipropi...	15238701	1
Caelo	Betamethasone dipropi...	15H28-B06	1
Caelo	Betamethasone dipropi...	15412906	10
Caelo	Betamethasone dipropi...	15412903	2
Caelo	Betamethasone dipropi...	18090606	5
Caelo	Betamethasone dipropi...	15412907	4
Caelo	Betamethasone dipropi...	17105907	3
Caelo	Betamethasone dipropi...	15163305	1
Caelo	Betamethasone dipropi...	15163304	1
Caelo	Betamethasone dipropi...	15412910	6
Caelo	Betamethasone dipropi...	17B16-B04-332220	1
Caelo	Betamethasone dipropi...	3170717	1
Caelo	Betamethasone dipropi...	15412913	1
Caelo	Betamethasone dipropi...	18090601	3
Caelo	Betamethasone dipropi...	17105903	1
Caelo	Betamethasone dipropi...	18090607	3
Caelo	Betamethasone dipropi...	17105901	9
Caelo	Betamethasone dipropi...	17105902	9
Caelo	Betamethasone dipropi...	15412916	1
Caelo	Betamethasone dipropi...	15412905	1
Caelo	Betamethasone dipropi...	15412909	1
Caelo	Betamethasone dipropi...	13023002	5
Caelo	Betamethasone dipropi...	150605BD	1
Caelo	Betamethasone dipropi...	18090602	2
Caelo	Betamethasone dipropi...	15412901	1
Caelo	Betamethasone dipropi...	15238701	21
Caelo	Betamethasone valerate	12003406	6
Caelo	Betamethasone valerate	15L21-B03-318663	1
Caelo	Betamethasone valerate	16270501	15
Caelo	Betamethasone valerate	15398502	1
Caelo	Betamethasone valerate	210112	1
Caelo	Betamethasone valerate	15398503	2
Caelo	Betamethasone valerate	15262612	1
Caelo	Betamethasone valerate	15398512	16
Caelo	Betamethasone valerate	15075811	1
Caelo	Betamethasone valerate	15262610	4
Caelo	Betamethasone valerate	12003405	8
Caelo	Betamethasone valerate	15398515	21
Caelo	Betamethasone valerate	BV018/0715	2
Caelo	Betamethasone valerate	160216BV	1
Caelo	Betamethasone valerate	15075804	7
Caelo	Betamethasone valerate	15398506	24
Caelo	Betamethasone valerate	1516-B07-319228	1
Caelo	Betamethasone valerate	161125BV/15398512	1
Caelo	Betamethasone valerate	1538506	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Betamethasone valerate	13091212	8
Caelo	Betamethasone valerate	15075805	6
Caelo	Betamethasone valerate	16720513	1
Caelo	Betamethasone valerate	12229810	3
Caelo	Betamethasone valerate	16270508	2
Caelo	Betamethasone valerate	16270513	15
Caelo	Betamethasone valerate	12003402	1
Caelo	Betamethasone valerate	12229812	9
Caelo	Betamethasone valerate	16270509	9
Caelo	Betamethasone valerate	16270511	16
Caelo	Betamethasone valerate	16270504	1
Caelo	Betamethasone valerate	13091201	4
Caelo	Betamethasone valerate	15398513	3
Caelo	Betamethasone valerate	16270503	13
Caelo	Betamethasone valerate	15398514	3
Caelo	Betamethasone valerate	16E23-B03	7
Caelo	Betamethasone valerate	14188605	1
Caelo	Betamethasone valerate	15398504	3
Caelo	Betamethasone valerate	1627051	1
Caelo	Betamethasone valerate	12229808	13
Caelo	Betamethasone valerate	B03041205	1
Caelo	Betamethasone valerate	15F16-B04-315327	1
Caelo	Betamethasone valerate	1610010-01	46
Caelo	Betamethasone valerate	13091206	1
Caelo	Betamethasone valerate	13012015C	1
Caelo	Betamethasone valerate	12229806	5
Caelo	Betamethasone valerate	1526213	1
Caelo	Betamethasone valerate	15075810	8
Caelo	Betamethasone valerate	12g09-n07	1
Caelo	Betamethasone valerate	15075802	20
Caelo	Betamethasone valerate	14188610	6
Caelo	Betamethasone valerate	12229809	2
Caelo	Betamethasone valerate	12229805	9
Caelo	Betamethasone valerate	14188617	2
Caelo	Betamethasone valerate	13091204	4
Caelo	Betamethasone valerate	11131909	1
Caelo	Betamethasone valerate	15075803	9
Caelo	Betamethasone valerate	10296701	1
Caelo	Betamethasone valerate	1488610	1
Caelo	Betamethasone valerate	15262603	11
Caelo	Betamethasone valerate	2x12229813	1
Caelo	Betamethasone valerate	10296705	2
Caelo	Betamethasone valerate	13091203	11
Caelo	Betamethasone valerate	12003403	9
Caelo	Betamethasone valerate	11131917	3
Caelo	Betamethasone valerate	1108045-03	10
Caelo	Betamethasone valerate	12229802	3
Caelo	Betamethasone valerate	1131914	1
Caelo	Betamethasone valerate	13091207	8
Caelo	Betamethasone valerate	14188619	2
Caelo	Betamethasone valerate	15262613	1
Caelo	Betamethasone valerate	18121206	1
Caelo	Betamethasone valerate	12229804	3
Caelo	Betamethasone valerate	12003410	1
Caelo	Betamethasone valerate	15262606	4
Caelo	Betamethasone valerate	15398501	16
Caelo	Betamethasone valerate	15262609	4
Caelo	Betamethasone valerate	11131905	1
Caelo	Betamethasone valerate	123345	1
Caelo	Betamethasone valerate	15/16-B07-319228	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Betamethasone valerate	12003407	1
Caelo	Betamethasone valerate	15262607	7
Caelo	Betamethasone valerate	13091209	6
Caelo	Betamethasone valerate	13K05_B08_292652	1
Caelo	Betamethasone valerate	1526604	1
Caelo	Betamethasone valerate	14188609	5
Caelo	Betamethasone valerate	150758036	1
Caelo	Betamethasone valerate	15075809	3
Caelo	Betamethasone valerate	15262604	9
Caelo	Betamethasone valerate	15262602	2
Caelo	Betamethasone valerate	12003411	4
Caelo	Betamethasone valerate	1113	1
Caelo	Betamethasone valerate	16270502	26
Caelo	Betamethasone valerate	12229811	6
Caelo	Betamethasone valerate	17271601	3
Caelo	Betamethasone valerate	1707815	2
Caelo	Betamethasone valerate	17271604	7
Caelo	Betamethasone valerate	17271603	14
Caelo	Betamethasone valerate	17078715	10
Caelo	Betamethasone valerate	17078720	12
Caelo	Betamethasone valerate	1707814	1
Caelo	Betamethasone valerate	17271602	1
Caelo	Betamethasone valerate	13225505	8
Caelo	Betamethasone valerate	17078719	13
Caelo	Betamethasone valerate	13225504	6
Caelo	Betamethasone valerate	13225508	10
Caelo	Betamethasone valerate	17078714	13
Caelo	Betamethasone valerate	16H22-B06-331273	1
Caelo	Betamethasone valerate	13225502	2
Caelo	Betamethasone valerate	13225515	7
Caelo	Betamethasone valerate	18071403	1
Caelo	Betamethasone valerate	13225507	2
Caelo	Betamethasone valerate	17329003	7
Caelo	Betamethasone valerate	17K16-B0/8-347310	1
Caelo	Betamethasone valerate	18199908	2
Caelo	Betamethasone valerate	18199916	2
Caelo	Betamethasone valerate	31(17329018)	1
Caelo	Betamethasone valerate	18199907	5
Caelo	Betamethasone valerate	1729002	1
Caelo	Betamethasone valerate	18199902	10
Caelo	Betamethasone valerate	18199906	2
Caelo	Betamethasone valerate	17329002	10
Caelo	Betamethasone valerate	17329011	2
Caelo	Betamethasone valerate	13091205	13
Caelo	Betamethasone valerate	17329019	1
Caelo	Betamethasone valerate	17329004	1
Caelo	Betamethasone valerate	17329017	2
Caelo	Betamethasone valerate	13225503	5
Caelo	Betamethasone valerate	17329015	5
Caelo	Betamethasone valerate	17329008	4
Caelo	Betamethasone valerate	17329018	23
Caelo	Betamethasone valerate	17C29-B04	4
Caelo	Betamethasone valerate	17329001	5
Caelo	Betamethasone valerate	12229813	2
Caelo	Betamethasone valerate	17/251/Bv	1
Caelo	Betamethasone valerate	17078103	1
Caelo	Betamethasone valerate	1708701	1
Caelo	Betamethasone valerate	16270514	3
Caelo	Betamethasone valerate	17078704	3
Caelo	Betamethasone valerate	17078705	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Betamethasone valerate	17078703	15
Caelo	Betamethasone valerate	13225516	10
Caelo	Betamethasone valerate	129711-5471	1
Caelo	Betamethasone valerate	170401901	6
Caelo	Betamethasone valerate	17078702	1
Caelo	Betamethasone valerate	13225607	1
Caelo	Betamethasone valerate	17078712	9
Caelo	Betamethasone valerate	14188602	3
Caelo	Betamethasone valerate	17078711	3
Caelo	Betamethasone valerate	13091213	7
Caelo	Betamethasone valerate	17078707	1
Caelo	Betamethasone valerate	17078713	12
Caelo	Betamethasone valerate	14188604	6
Caelo	Betamethasone valerate	14188603	6
Caelo	Betamethasone valerate	13225509	2
Caelo	Betamethasone, micron. ...	11/14-N02	2
Caelo	Betamethasone-17-vale. ...	12A05-N08	3
Caelo	Betamethasone-17-vale. ...	12a05-no8	1
Caelo	Betony herb, cut	13367802	1
Caelo	Betony herb, cut	13048906	1
Caelo	Betony herb, cut	134048903	1
Caelo	Betony herb, cut	12007201	1
Caelo	Betony herb, cut	12007208	1
Caelo	Bifonazole	14k10-B01-312067	1
Caelo	Bifonazole	18009702	3
Caelo	Bifonazole	15152701	3
Caelo	Bifonazole	1703041-01	4
Caelo	Bilberry fruit, dried	13130501	1
Caelo	Bilberry fruit, dried	18101307	1
Caelo	Bilberry fruit, dried	13130504	1
Caelo	Bilberry fruit, dried	26061301	1
Caelo	Bilberry fruit, dried	10265601	1
Caelo	Bilberry fruit, dried	10273708	1
Caelo	Bilberry fruit, dried	10273711	1
Caelo	Bilberry fruit, dried	12183701	2
Caelo	Bilberry fruit, dried	10273709	1
Caelo	Bilberry fruit, dried	10265603	1
Caelo	Bilberry fruit, dried	12183704	4
Caelo	Bilberry leaf, cut	12140204	4
Caelo	Bilberry leaf, cut	23111201	1
Caelo	Bilberry leaf, cut	9367603	2
Caelo	Bilberry leaf, cut	13009904	3
Caelo	Bilberry leaf, cut	3091302	1
Caelo	Bilberry leaf, cut	12140206	2
Caelo	Bilberry leaf, cut	1300907	1
Caelo	Bilberry leaf, cut	13009907	2
Caelo	Bilberry leaf, cut	13009901	3
Caelo	Bilberry leaf, cut	12140201	2
Caelo	Bilberry leaf, cut	13009903	4
Caelo	Biotin	12A23-N03	2
Caelo	Birch leaf, cut	11189508	1
Caelo	Birch leaf, cut	11189510	3
Caelo	Birch leaf, cut	11189502	1
Caelo	Birch leaf, cut	12193004	5
Caelo	Birch leaf, cut	107063	1
Caelo	Birch leaf, cut	12193105	3
Caelo	Birch leaf, cut	12193104	1
Caelo	Birch leaf, cut	11189507	1
Caelo	Birch leaf, cut	12193109	15
Caelo	Birch leaf, cut	13108906	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Birch leaf, cut	13108917	1
Caelo	Birch leaf, cut	14360702	1
Caelo	Birch leaf, cut	13337302	1
Caelo	Birch leaf, cut	13108916	1
Caelo	Birch leaf, cut	13337303-3561412013	1
Caelo	Birch leaf, cut	13337303	1
Caelo	Birch leaf, cut	13108905	8
Caelo	Birch leaf, cut	12193108	1
Caelo	Birch leaf, cut	IntCh:311092	1
Caelo	Bismuth gallath, basic	15200904	3
Caelo	Bismuth gallath, basic	13179111	2
Caelo	Bismuth gallath, basic	16038312	6
Caelo	Bismuth gallath, basic	17-E16-B02-336992	1
Caelo	Bismuth gallath, basic	1520910	1
Caelo	Bismuth gallath, basic	2017002570	1
Caelo	Bismuth gallath, basic	17004102	7
Caelo	Bismuth gallath, basic	17004106	7
Caelo	Bismuth gallath, basic	62183516	1
Caelo	Bismuth gallath, basic	81297468	1
Caelo	Bismuth gallath, basic	13I18-B13-294108	1
Caelo	Bismuth gallath, basic	12172412	3
Caelo	Bismuth gallath, basic	15080602	2
Caelo	Bismuth gallath, basic	14224103	3
Caelo	Bismuth gallath, basic	1520913	2
Caelo	Bismuth gallath, basic	15J26_B02-321313	1
Caelo	Bismuth gallath, basic	9K26-N10	2
Caelo	Bismuth gallath, basic	12172402	1
Caelo	Bismuth gallath, basic	14224103	2
Caelo	Bismuth gallath, basic	15J26-B02-320479	1
Caelo	Bismuth gallath, basic	15200907	1
Caelo	Bismuth gallath, basic	15B02B04304413	1
Caelo	Bismuth gallath, basic	11127104	1
Caelo	Bismuth gallath, basic	23121410	1
Caelo	Bismuth gallath, basic	18A29-B06	1
Caelo	Bismuth gallath, basic	16038308	7
Caelo	Bismuth gallath, basic	17259503	1
Caelo	Bismuth gallath, basic	15080603	1
Caelo	Bismuth gallath, basic	17259509	3
Caelo	Bismuth gallath, basic	10135109	1
Caelo	Bismuth gallath, basic	15200905	1
Caelo	Bismuth gallath, basic	15200913	4
Caelo	Bismuth gallath, basic	18A29-B06-346305	3
Caelo	Bismuth gallath, basic	18040801	3
Caelo	Bismuth gallath, basic	12172409	1
Caelo	Bismuth gallath, basic	1725909	1
Caelo	Bismuth gallath, basic	17004112	6
Caelo	Bismuth gallath, basic	Nachtestat04,06,2016	1
Caelo	Bismuth gallath, basic	81297468	1
Caelo	Bismuth gallath, basic	17004109	3
Caelo	Bismuth gallath, basic	15B02-B04-309481	5
Caelo	Bismuth gallath, basic	13179105	1
Caelo	Bismuth gallath, basic	15080602	2
Caelo	Bismuth gallath, basic	15200911	6
Caelo	Bismuth gallath, basic	15080606	9
Caelo	Bismuth gallath, basic	B-150513-1	1
Caelo	Bismuth gallath, basic	92057409	1
Caelo	Bismuth gallath, basic	12166901	2
Caelo	Bismuth gallath, basic	13179111	2
Caelo	Bismuth gallath, basic	12172409	5
Caelo	Bismuth gallath, basic	12I13-N11	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Bismuth gallath, basic	12C13-N03	5
Caelo	Bismuth gallath, basic	1127203	1
Caelo	Bismuth gallath, basic	1217240	1
Caelo	Bismuth gallath, basic	8129768	1
Caelo	Bismuth gallath, basic	9205703	2
Caelo	Bismuth gallath, basic	12172402	4
Caelo	Bismuth gallath, basic	13/18-B13-292057	2
Caelo	Bismuth gallath, basic	13179101	1
Caelo	Bismuth gallath, basic	5208I-01185	1
Caelo	Bismuth gallath, basic	12c13-n03	1
Caelo	Bismuth gallath, basic	12172401	2
Caelo	Bismuth gallath, basic	12c13-no3	1
Caelo	Bismuth gallath, basic	9205701	1
Caelo	Bismuth gallath, basic	81297468	3
Caelo	Bismuth gallath, basic	2414E-01185	1
Caelo	Bismuth gallath, basic	10D23-N05	1
Caelo	Bitter fennel	1S443206	1
Caelo	Bitter fennel	15443212	3
Caelo	Bitter fennel	12094710	4
Caelo	Bitter fennel	13200402	2
Caelo	Bitter fennel	13148306	1
Caelo	Bitter fennel	15443202	3
Caelo	Bitter fennel	12094721	3
Caelo	Bitter fennel	15006508	2
Caelo	Bitter fennel	13148305	2
Caelo	Bitter fennel	13148304-3500411070	1
Caelo	Bitter fennel	13371610	1
Caelo	Bitter fennel	12094705	3
Caelo	Bitter fennel	6091201	1
Caelo	Bitter fennel	12094712	4
Caelo	Bitter fennel	12094704	3
Caelo	Bitter fennel	5061302	1
Caelo	Bitter fennel	12094719	4
Caelo	Bitter fennel	13398502	1
Caelo	Bitter fennel	15006506	1
Caelo	Bitter fennel	15443203	2
Caelo	Bitter fennel	12035102	1
Caelo	Bitter fennel	402311545	1
Caelo	Bitter fennel	11016330	1
Caelo	Bitter fennel	14376401	1
Caelo	Bitter fennel	12094725	5
Caelo	Bitter fennel	12094707	4
Caelo	Bitter fennel	15443204	7
Caelo	Bitter fennel	14255802	3
Caelo	Bitter fennel	15443208	2
Caelo	Bitter fennel	15443205	2
Caelo	Bitter fennel	16005902	1
Caelo	Bitter fennel	12167102	4
Caelo	Bitter fennel	16061202	1
Caelo	Bitter fennel	5121401	1
Caelo	Bitter fennel	15443209	3
Caelo	Bitter fennel	13371605	1
Caelo	Bitter fennel	12035102	2
Caelo	Bitter fennel	1302081	1
Caelo	Bitter fennel	12094716	6
Caelo	Bitter fennel	13148307	1
Caelo	Bitter fennel	7031306	1
Caelo	Bitter fennel	14376402	1
Caelo	Bitter fennel	14376408;21S08	1
Caelo	Bitter fennel	2150317	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Bitter fennel	12094724	3
Caelo	Bitter fennel	17020603	1
Caelo	Bitter fennel	12167101	1
Caelo	Bitter fennel	12362710	1
Caelo	Bitter fennel	15006501	1
Caelo	Bitter fennel	16298102	1
Caelo	Bitter fennel	16298101	1
Caelo	Bitter fennel	1320041	1
Caelo	Bitter fennel	13148303	2
Caelo	Bitter fennel	13020801	4
Caelo	Bitter fennel	14111301	1
Caelo	Bitter fennel	12094703	2
Caelo	Bitter fennel	13148304	18
Caelo	Bitter fennel	12094716-3408310191	1
Caelo	Bitter fennel	1209474	1
Caelo	Bitter fennel	12094715	4
Caelo	Bitter fennel	12345678	1
Caelo	Bitter fennel	15443213	1
Caelo	Bitter fennel	11185501	1
Caelo	Bitter fennel	12094726	2
Caelo	Bitter fennel	4061302	1
Caelo	Bitter fennel	14376406	1
Caelo	Bitter fennel	12094717	1
Caelo	Bitter fennel	13398502	1
Caelo	Bitter fennel	12094720	5
Caelo	Bitter fennel	15007601	1
Caelo	Bitter fennel	15141004	1
Caelo	Bitter fennel	13371610	2
Caelo	Bitter fennel	27121201	1
Caelo	Bitter fennel	154443204	1
Caelo	Bitter fennel	15443215	4
Caelo	Bitter fennel	14376405	1
Caelo	Bitter fennel	9121305	1
Caelo	Bitter fennel	15443216	1
Caelo	Bitter fennel	13200401	4
Caelo	Bitter fennel	15006507	1
Caelo	Bitter fennel	1250117/15443216	1
Caelo	Bitter fennel	13020802	6
Caelo	Bitter orange flower ...	11049202	1
Caelo	Bitter orange flower ...	12266703	2
Caelo	Bitter orange flower ...	12266605	1
Caelo	Bitter orange flower ...	12266608	6
Caelo	Bitter orange flower ...	11049212	1
Caelo	Bitter orange flower ...	12266705	1
Caelo	Bitter orange flower ...	12266609	1
Caelo	Bitter orange flower ...	13352701-3518411022	1
Caelo	Bitter orange flower ...	12266614	1
Caelo	Bitter orange flower ...	2712E-01445	1
Caelo	Bitter orange flower ...	12266605	1
Caelo	Bitter orange flower ...	12266707	1
Caelo	Bitter orange flower ...	12266611	1
Caelo	Bitter orange flower ...	13067903	2
Caelo	Bitter orange flower ...	12266708	1
Caelo	Bitter orange flower ...	12266611	1
Caelo	Bitter orange flower ...	12266608	2
Caelo	Bitter orange flower ...	6051301	1
Caelo	Bitter orange flower ...	13352701	1
Caelo	Bitter orange peel, cut	12067402	1
Caelo	Bitter orange peel, cut	12067405	1
Caelo	Bitter orange peel, cut	27011403	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Bitter orange peel, cut	12239203	3
Caelo	Bitter orange peel, cut	4312M-01531	1
Caelo	Bitter orange peel, cut	21101302	1
Caelo	Bitter orange peel, cut	17111201	1
Caelo	Bitter orange peel, cut	13049001	1
Caelo	Black cohosh, cut	10281104	1
Caelo	Black cohosh, cut	13237903	1
Caelo	Black cohosh, cut	13237901	2
Caelo	Black cohosh, cut	10281118	2
Caelo	Black cohosh, cut	13142101	1
Caelo	Black cohosh, cut	10281112	1
Caelo	Black cumin seed	13122013b	1
Caelo	Black cumin seed	16115708	1
Caelo	Black cumin seed	13326406	1
Caelo	Black cumin seed	16115703	1
Caelo	Black cumin seed	15300401	1
Caelo	Black cumin seed	14359905	1
Caelo	Black cumin seed	10121915	1
Caelo	Black cumin seed	11233507	1
Caelo	Black cumin seed	12358906	2
Caelo	Black cumin seed	12358902	1
Caelo	Black cumin seed	11233503	1
Caelo	Black cumin seed	13326405	2
Caelo	Black cumin seed	11233502	2
Caelo	Black cumin seed	14359903	3
Caelo	Black cumin seed	12358904	1
Caelo	Black cumin seed	12358901	2
Caelo	Black cumin seed	12358903	2
Caelo	Black cumin seed	10121916	1
Caelo	Black mustard seed	13038904	1
Caelo	Black mustard seed	12189307	1
Caelo	Black mustard seed	14057702	1
Caelo	Black mustard seed, p...	12009703	1
Caelo	Black mustard seed, p...	11276701	1
Caelo	Black mustard seed, p...	12009706	1
Caelo	Black mustard seed, p...	12009702	4
Caelo	Black mustard seed, p...	1633503	1
Caelo	Black mustard seed, p...	12208604	3
Caelo	Black mustard seed, p...	14162808	1
Caelo	Black mustard seed, p...	14293710	2
Caelo	Black mustard seed, p...	14293706	1
Caelo	Black mustard seed, p...	14293713	1
Caelo	Black mustard seed, p...	2015Q-02244	1
Caelo	Black mustard seed, p...	15131104	1
Caelo	Black mustard seed, p...	15131102	3
Caelo	Black mustard seed, p...	15131105	1
Caelo	Black mustard seed, p...	15131101	1
Caelo	Black mustard seed, p...	15214703	4
Caelo	Black mustard seed, p...	15214706	3
Caelo	Black mustard seed, p...	15214703	1
Caelo	Black mustard seed, p...	15214708	4
Caelo	Black mustard seed, p...	14162802	1
Caelo	Black mustard seed, p...	15214713	2
Caelo	Black mustard seed, p...	15214714	1
Caelo	Black mustard seed, p...	16051802	2
Caelo	Black mustard seed, p...	16051807	1
Caelo	Black mustard seed, p...	16051807	1
Caelo	Black mustard seed, p...	16051810	3
Caelo	Black mustard seed, p...	16272704	2
Caelo	Black mustard seed, p...	16272703	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Black mustard seed, p...	16272707	3
Caelo	Black mustard seed, p...	17276404	1
Caelo	Black mustard seed, p...	17196303	4
Caelo	Black mustard seed, p...	17196301	1
Caelo	Black mustard seed, p...	17196306	1
Caelo	Black mustard seed, p...	16333501	2
Caelo	Black mustard seed, p...	12208602	5
Caelo	Black mustard seed, p...	14293702	6
Caelo	Black mustard seed, p...	15214707	1
Caelo	Black mustard seed, p...	15214707	1
Caelo	Black mustard seed, p...	12350207	5
Caelo	Black mustard seed, p...	12350202-3369319012	1
Caelo	Black mustard seed, p...	13097702	4
Caelo	Black mustard seed, p...	12350204	2
Caelo	Black mustard seed, p...	12350210	1
Caelo	Black mustard seed, p...	12350211	3
Caelo	Black mustard seed, p...	12350208	7
Caelo	Black mustard seed, p...	13097706	2
Caelo	Black mustard seed, p...	16333503	2
Caelo	Black mustard seed, p...	13097705	6
Caelo	Black mustard seed, p...	12350209	1
Caelo	Black mustard seed, p...	12350205	2
Caelo	Black mustard seed, p...	14293701	2
Caelo	Black mustard seed, p...	13097701	11
Caelo	Black mustard seed, p...	12350202	16
Caelo	Black mustard seed, p...	13358302	13
Caelo	Black mustard seed, p...	12208601	3
Caelo	Black mustard seed, p...	25011401	1
Caelo	Black mustard seed, p...	13358304	4
Caelo	Black mustard seed, p...	23061401	1
Caelo	Black mustard seed, p...	14162810	2
Caelo	Black mustard seed, p...	14162810	3
Caelo	Black mustard seed, p...	14162807	2
Caelo	Black mustard seed, p...	14245601	1
Caelo	Black mustard seed, p...	13097704	4
Caelo	Black pepper	12192309	1
Caelo	Black pepper	12192305	1
Caelo	Black pepper	13221110	1
Caelo	Black pepper	12192306	1
Caelo	Black pepper	15395310	1
Caelo	Black pepper	12192304	1
Caelo	Black pepper	11193911	1
Caelo	Blackberry leaf, cut	11030704	4
Caelo	Blackberry leaf, cut	10136403	1
Caelo	Blackberry leaf, cut	10136408	6
Caelo	Blackberry leaf, cut	28061202	1
Caelo	Blackberry leaf, cut	13091206	1
Caelo	Blackberry leaf, cut	11030701	2
Caelo	Blackberry leaf, cut	11030707	5
Caelo	Blackberry leaf, cut	3011E-01512	1
Caelo	Blackberry leaf, cut	11030706	1
Caelo	Blackberry leaf, cut	3091303	1
Caelo	Blackberry leaf, cut	2271c120508	1
Caelo	Blackberry leaf, cut	15039301	1
Caelo	Blackberry leaf, cut	14111318	1
Caelo	Blackberry leaf, cut	14019010	1
Caelo	Blackberry leaf, cut	13127901	1
Caelo	Blackberry leaf, cut	13101610	1
Caelo	Blackcurrant leaf, cut	12168804	1
Caelo	Blackcurrant leaf, cut	9373913	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Blackcurrant leaf, cut	12115902	1
Caelo	Blackcurrant leaf, cut	12168810	3
Caelo	Blackcurrant leaf, cut	80558228	1
Caelo	Blackcurrants	12272804	3
Caelo	Blackthorn flower	12178105	1
Caelo	Blackthorn flower	12178103	2
Caelo	Blackthorn flower	11206906	1
Caelo	Blackthorn flower	12178101	1
Caelo	Blackthorn flower	12178104	4
Caelo	Blackthorn flower	12178102	9
Caelo	Blackthorn flower, ru...	1150N101109	2
Caelo	Blessed thistle, cut	13025902	1
Caelo	Blessed thistle, cut	13025901	5
Caelo	Blessed thistle, cut	10325404	1
Caelo	Blessed thistle, cut	27111205	1
Caelo	Blessed thistle, cut	13025903	3
Caelo	Blessed thistle, cut	4743554	1
Caelo	Blessed thistle, cut	3091301	1
Caelo	Blessed thistle, cut	12038202	1
Caelo	Blessed thistle, cut	9355208	1
Caelo	Blessed thistle, cut	10325402	5
Caelo	Blessed thistle, cut	12038203	2
Caelo	Blessed thistle, cut	13389701	1
Caelo	Blue bottle flower wi...	12232501	2
Caelo	Blue bottle flower wi...	15206601	1
Caelo	Blue bottle flower wi...	14323302	1
Caelo	Blue bottle flower wi...	12060205	4
Caelo	Blue bottle flower wi...	11059211	2
Caelo	Blue bottle flower wi...	12060203	4
Caelo	Blue bottle flower wi...	12060202	2
Caelo	Blue bottle flower wi...	11059206	2
Caelo	Blue bottle flower wi...	11059209	2
Caelo	Bogbean leaf, cut	12226101	1
Caelo	Bogbean leaf, cut	1226101	1
Caelo	Bogbean leaf, cut	12226104	1
Caelo	Bogbean leaf, cut	13123901	2
Caelo	Bogbean leaf, cut	13149001	1
Caelo	Bogbean leaf, cut	13123902	2
Caelo	Boldo leaf, cut	655456	1
Caelo	Boldo leaf, cut	12064303	1
Caelo	Boldo leaf, cut	12064305	1
Caelo	Boldo leaf, cut	11016511	1
Caelo	Boldo leaf, cut	11016509	1
Caelo	Boldo leaf, cut	12227005	1
Caelo	Boldo leaf, cut	12227004	1
Caelo	Boldo leaf, cut	12064306	1
Caelo	Boric acid	15425101	5
Caelo	Boric acid	80303AB	1
Caelo	Boric acid	15158903	1
Caelo	Boric acid	9324021	1
Caelo	Boric acid	9324043	1
Caelo	Boric acid	14133119	1
Caelo	Boric acid	15095103	2
Caelo	Boric acid	15425108	2
Caelo	Boric acid	12238901	1
Caelo	Boric acid	14095202	1
Caelo	Boric acid	15158901	1
Caelo	Boric acid	12260502	1
Caelo	Boric acid	10227913	1
Caelo	Boric acid	15425112	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Boric acid	14133113	2
Caelo	Boric acid	15242909	1
Caelo	Boric acid	17327907	2
Caelo	Boric acid	15425114	1
Caelo	Boric acid	15425118	2
Caelo	Boric acid	15242906	1
Caelo	Boric acid	5425114	1
Caelo	Boric acid	9324036	1
Caelo	Boric acid	11D13-N08	1
Caelo	Bucho leaf	13056602	2
Caelo	Bucho leaf	13056605	1
Caelo	Bucho leaf	13056601	2
Caelo	Buckwheat herb, cut	5111302	1
Caelo	Buckwheat herb, cut	30041301	1
Caelo	Buckwheat herb, cut	13193203	1
Caelo	Buckwheat herb, cut	11212810	1
Caelo	Buckwheat herb, cut	11212612	1
Caelo	Buckwheat herb, cut	11212807	2
Caelo	Buckwheat herb, cut	11212817	6
Caelo	Buckwheat herb, cut	11212808	2
Caelo	Buckwheat herb, cut	11212814	2
Caelo	Buckwheat herb, cut	11212809	1
Caelo	Buckwheat herb, cut	13051303	1
Caelo	Buckwheat herb, cut	2051303	1
Caelo	Buckwheat herb, cut	11212812	1
Caelo	Budesonide, micronized	14071701	3
Caelo	Budesonide, micronized	12295108	1
Caelo	Budesonide, micronized	12295106	2
Caelo	Budesonide, micronized	12295105	1
Caelo	Budesonide, micronized	11175015	1
Caelo	Budesonide, micronized	12295109	2
Caelo	Budesonide, micronized	12295103	1
Caelo	Budesonide, micronized	17078801	8
Caelo	Budesonide, micronized	11351703	1
Caelo	Budesonide, micronized	11351707	3
Caelo	Budesonide, micronized	21100M1	1
Caelo	Budesonide, micronized	11351701	1
Caelo	Budesonide, micronized	11351706	5
Caelo	Budesonide, micronized	16143106	1
Caelo	Budesonide, micronized	15344012	4
Caelo	Budesonide, micronized	12295110	4
Caelo	Budesonide, micronized	14071717	3
Caelo	Budesonide, micronized	15344008	6
Caelo	Budesonide, micronized	15344013	1
Caelo	Budesonide, micronized	15344010	3
Caelo	Budesonide, micronized	15344004	1
Caelo	Budesonide, micronized	17078805	6
Caelo	Budesonide, micronized	15044002	1
Caelo	Budesonide, micronized	15344014	1
Caelo	Budesonide, micronized	16143102	2
Caelo	Budesonide, micronized	15344002	1
Caelo	Budesonide, micronized	15173802	2
Caelo	Budesonide, micronized	15173807	3
Caelo	Budesonide, micronized	12295111	2
Caelo	Budesonide, micronized	15173811	6
Caelo	Budesonide, micronized	15173810	1
Caelo	Budesonide, micronized	16143112	2
Caelo	Budesonide, micronized	15344003	6
Caelo	Budesonide, micronized	16143103	3
Caelo	Budesonide, micronized	16143101	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Budesonide, micronized	15173803	3
Caelo	Budesonide, micronized	15173806	2
Caelo	Budesonide, micronized	774865	1
Caelo	Budesonide, micronized	1471707	1
Caelo	Budesonide, micronized	14071714	4
Caelo	Budesonide, micronized	16143108	7
Caelo	Budesonide, micronized	15344005	2
Caelo	Budesonide, micronized	15344006	1
Caelo	Budesonide, micronized	11351713	1
Caelo	Budesonide, micronized	17342902	1
Caelo	Budesonide, micronized	18172610	1
Caelo	Budesonide, micronized	18172611	1
Caelo	Budesonide, micronized	18172604	1
Caelo	Budesonide, micronized	18074909	1
Caelo	Budesonide, micronized	18172603	4
Caelo	Budesonide, micronized	18074907	2
Caelo	Budesonide, micronized	18074905	1
Caelo	Budesonide, micronized	18074902	2
Caelo	Budesonide, micronized	11175014	2
Caelo	Budesonide, micronized	18074903	1
Caelo	Budesonide, micronized	18074901	3
Caelo	Budesonide, micronized	17078808	5
Caelo	Budesonide, micronized	17078803	2
Caelo	Budesonide, micronized	17078809	3
Caelo	Budesonide, micronized	170078801	1
Caelo	Budesonide, micronized	16143105	2
Caelo	Budesonide, micronized	17078804	2
Caelo	Budesonide, micronized	16143107	3
Caelo	Burdock root, cut	12295808	2
Caelo	Burdock root, cut	12295801	2
Caelo	Burdock root, cut	107063	1
Caelo	Burdock root, cut	12295807	1
Caelo	Burdock root, cut	13122802	1
Caelo	Burdock root, cut	11100801	2
Caelo	Burdock root, cut	13122801	2
Caelo	Burdock root, cut	10073213	1
Caelo	Burdock root, cut	12295805	3
Caelo	Burdock root, cut	1100801	1
Caelo	Burdock root, cut	12295802	4
Caelo	Burnet root	8111203	1
Caelo	Burnet root	12027403	1
Caelo	Burnet root	12027405	2
Caelo	Burnet root	12027413	1
Caelo	Burnet root	12027409	1
Caelo	Burnet root	12027410	3
Caelo	Burnet root	12027411	3
Caelo	Butylated hydroxytolu...	14300503	2
Caelo	Butylated hydroxytolu...	15172105	10
Caelo	Butylated hydroxytolu...	14300501	9
Caelo	Butylated hydroxytolu...	13354604	3
Caelo	Butylated hydroxytolu...	14300502	5
Caelo	Butylated hydroxytolu...	12280206	1
Caelo	Butylated hydroxytolu...	15172103	20
Caelo	Butylated hydroxytolu...	143004502	1
Caelo	Butylated hydroxytolu...	14300506	1
Caelo	Butylated hydroxytolu...	14300505	6
Caelo	Butylated hydroxytolu...	15094501	4
Caelo	Butylated hydroxytolu...	15094501	4
Caelo	Butylated hydroxytolu...	1400501	1
Caelo	Butylated hydroxytolu...	15172101	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Butylated hydroxytolu...	18029104	4
Caelo	Butylated hydroxytolu...	18029107	7
Caelo	Butylated hydroxytolu...	2113U02851	1
Caelo	Butylated hydroxytolu...	14300505	1
Caelo	Butylated hydroxytolu...	18029105	13
Caelo	Butylated hydroxytolu...	19074902	1
Caelo	Butylated hydroxytolu...	12280206	3
Caelo	Butylated hydroxytolu...	122080206	1
Caelo	Butylated hydroxytolu...	12207005	3
Caelo	Butylated hydroxytolu...	163163055402	1
Caelo	Butylated hydroxytolu...	16316305	11
Caelo	Butylated hydroxytolu...	18229102	1
Caelo	Butylated hydroxytolu...	16316304-20180619Ma5	1
Caelo	Butylated hydroxytolu...	20161112-1	1
Caelo	Butylated hydroxytolu...	1631301	1
Caelo	Butylated hydroxytolu...	1220705	1
Caelo	Butylated hydroxytolu...	18029102	11
Caelo	Butylated hydroxytolu...	15426807	12
Caelo	Butylated hydroxytolu...	12280201	1
Caelo	Butylated hydroxytolu...	12280203	5
Caelo	Butylated hydroxytolu...	12382302	5
Caelo	Butylated hydroxytolu...	123823	1
Caelo	Butylated hydroxytolu...	12382304	8
Caelo	Butylated hydroxytolu...	12382309	8
Caelo	Butylated hydroxytolu...	18029108	1
Caelo	Butylated hydroxytolu...	15172101	8
Caelo	Butylated hydroxytolu...	16316307	16
Caelo	Butylated hydroxytolu...	16316304	12
Caelo	Butylated hydroxytolu...	11(18029107)	1
Caelo	Butylated hydroxytolu...	231115BHT	1
Caelo	Butylated hydroxytolu...	11131006	1
Caelo	Butylated hydroxytolu...	13354611	3
Caelo	Butylated hydroxytolu...	13354604	6
Caelo	Butylated hydroxytolu...	12382307	1
Caelo	Butylated hydroxytolu...	14300501	1
Caelo	Butylated hydroxytolu...	12382312	1
Caelo	Butylated hydroxytolu...	18048701	12
Caelo	Butylated hydroxytolu...	18029105	1
Caelo	Butylated hydroxytolu...	1426801	1
Caelo	Butylated hydroxytolu...	716M-02851	1
Caelo	Butylated hydroxytolu...	16316301	10
Caelo	Butylated hydroxytolu...	7031403	1
Caelo	Butylated hydroxytolu...	12158204	5
Caelo	Butylated hydroxytolu...	13354602	4
Caelo	Butylated hydroxytolu...	2170717	1
Caelo	Caffeine	18011906	3
Caelo	Caffeine	11076807	1
Caelo	Caffeine	15372012	1
Caelo	Caffeine	18011902	2
Caelo	Caffeine	18011901	4
Caelo	Caffeine	15082311	1
Caelo	Caffeine	15375006	2
Caelo	Caffeine	15375007	2
Caelo	Caffeine	11111111	1
Caelo	Caffeine	18077906	1
Caelo	Caffeine	17173704	4
Caelo	Caffeine	15375012	4
Caelo	Caffeine	17173707	1
Caelo	Caffeine	15082315	1
Caelo	Caffeine	15054709	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Caffeine	1211012-01	1
Caelo	Caffeine	12317909	7
Caelo	Caffeine	3212E-01281	1
Caelo	Caffeine	12093418	2
Caelo	Caffeine	12317906	1
Caelo	Caffeine	12317904	2
Caelo	Caffeine	12093410	2
Caelo	Caffeine	12317914	2
Caelo	Caffeine	13233105	2
Caelo	Caffeine	9301807	1
Caelo	Caffeine	13233105	2
Caelo	Caffeine	12317902	1
Caelo	Caffeine	18011911	2
Caelo	Caffeine	12093402	2
Caelo	Caffeine	13233111	3
Caelo	Caffeine	13233111	2
Caelo	Caffeine	15054704	2
Caelo	Caffeine	1503020-01	8
Caelo	Caffeine	11076807	1
Caelo	Caffeine	18011915	4
Caelo	Caffeine citrate	13D30-N02	9
Caelo	Caffeine citrate	11F21-N05	1
Caelo	Calcium ascorbate	15324708	1
Caelo	Calcium ascorbate	13105608	1
Caelo	Calcium ascorbate	15324702	1
Caelo	Calcium ascorbate	16311105	1
Caelo	Calcium ascorbate	14340907	1
Caelo	Calcium chloride dihy...	12341406	1
Caelo	Calcium chloride dihy...	8237314	1
Caelo	Calcium chloride dihy...	12341409	1
Caelo	Calcium chloride dihy...	12341404	2
Caelo	Calcium chloride dihy...	11230915	1
Caelo	Calcium chloride dihy...	11230921	2
Caelo	Calcium chloride dihy...	12B09-N11	1
Caelo	Calcium citrate tetra...	15177701	1
Caelo	Calcium citrate tetra...	11168806	3
Caelo	Calcium citrate tetra...	1509039-01	4
Caelo	Calcium citrate tetra...	15375702	6
Caelo	Calcium citrate tetra...	15375707	1
Caelo	Calcium citrate tetra...	15177703	7
Caelo	Calcium citrate tetra...	11349325Nachprüfung	1
Caelo	Calcium citrate tetra...	14133310	1
Caelo	Calcium citrate tetra...	11349318	3
Caelo	Calcium citrate tetra...	15177702	2
Caelo	Calcium citrate tetra...	11349320	1
Caelo	Calcium citrate tetra...		1
Caelo	Calcium citrate tetra...	12310402	2
Caelo	Calcium citrate tetra...	12310406	5
Caelo	Calcium citrate tetra...	71356168	1
Caelo	Calcium citrate tetra...	1168806	1
Caelo	Calcium citrate tetra...	19061315	1
Caelo	Calcium citrate tetra...	12310411	2
Caelo	Calcium citrate tetra...	1168813	1
Caelo	Calcium citrate tetra...	12310413	4
Caelo	Calcium citrate tetra...	12310407	2
Caelo	Calcium citrate tetra...	92457399	1
Caelo	Calcium citrate tetra...	12310417	3
Caelo	Calcium citrate tetra...	12310403	3
Caelo	Calcium citrate tetra...	12310418	2
Caelo	Calcium citrate tetra...	133912042x	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Calcium citrate tetra...	13391205	5
Caelo	Calcium citrate tetra...	13391203	4
Caelo	Calcium citrate tetra...	2114M1209	1
Caelo	Calcium citrate tetra...	3512U-01209	1
Caelo	Calcium citrate tetra...	14133303	1
Caelo	Calcium citrate tetra...	16021502	1
Caelo	Calcium citrate tetra...	11349325	2
Caelo	Calcium citrate tetra...	15375712	2
Caelo	Calcium citrate tetra...	15375703	1
Caelo	Calcium citrate tetra...	15375713	6
Caelo	Calcium citrate tetra...	19010002-01	1
Caelo	Calcium citrate tetra...	17176812	2
Caelo	Calcium citrate tetra...	11168813	1
Caelo	Calcium citrate tetra...	17176805	1
Caelo	Calcium citrate tetra...	17176806	1
Caelo	Calcium citrate tetra...	17176808	7
Caelo	Calcium citrate tetra...	1716803	1
Caelo	Calcium citrate tetra...	11349304	2
Caelo	Calcium citrate tetra...	17176803	5
Caelo	Calcium citrate tetra...	16056510	2
Caelo	Calcium citrate tetra...	1600565	1
Caelo	Calcium citrate tetra...	11349309	2
Caelo	Calcium citrate tetra...	11349301	1
Caelo	Calcium citrate tetra...	11349313	2
Caelo	Calcium citrate tetra...	92192429	1
Caelo	Calcium citrate tetra...	16056505	2
Caelo	Calcium citrate tetra...	16056504	4
Caelo	Calcium citrate tetra...	16056501	1
Caelo	Calcium citrate tetra...	11168803	1
Caelo	Calcium gluconate	15126301	1
Caelo	Calcium gluconate	18101003	2
Caelo	Calcium gluconate	4(18101004)	1
Caelo	Calcium gluconate	1606034	1
Caelo	Calcium gluconate	15435611	1
Caelo	Calcium gluconate	18101010	1
Caelo	Calcium gluconate	12335813	1
Caelo	Calcium gluconate	15435607	2
Caelo	Calcium gluconate	15324807	1
Caelo	Calcium gluconate	18101004	4
Caelo	Calcium gluconate	15435620	2
Caelo	Calcium gluconate	1000041-01	1
Caelo	Calcium gluconate	14085405	1
Caelo	Calcium gluconate	15435619	1
Caelo	Calcium gluconate	16603401	1
Caelo	Calcium gluconate	15435613	1
Caelo	Calcium gluconate	15435601	2
Caelo	Calcium gluconate	14230706	1
Caelo	Calcium gluconate	18106202	1
Caelo	Calcium gluconate	15435622	2
Caelo	Calcium gluconate	12335810	1
Caelo	Calcium gluconate	15324806	1
Caelo	Calcium gluconate	1543601	1
Caelo	Calcium glycerophosph...	17250703	2
Caelo	Calcium glycerophosph...	16231007	2
Caelo	Calcium glycerophosph...	16231001	3
Caelo	Calcium glycerophosph...	162310	1
Caelo	Calcium glycerophosph...	16231006	3
Caelo	Calcium glycerophosph...	14341013	1
Caelo	Calcium glycerophosph...	12012413	2
Caelo	Calcium glycerophosph...	16231008	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Calcium glycerophosph...	16231005	1
Caelo	Calcium glycerophosph...	16231010	3
Caelo	Calcium glycerophosph...	17359101	5
Caelo	Calcium glycerophosph...	12012406	1
Caelo	Calcium glycerophosph...	14341015	1
Caelo	Calcium glycerophosph...	14347101	2
Caelo	Calcium glycerophosph...	15B03-B03-305050	3
Caelo	Calcium glycerophosph...	14341010	1
Caelo	Calcium glycerophosph...	14341007	5
Caelo	Calcium glycerophosph...	19100703	1
Caelo	Calcium glycerophosph...	13007914	1
Caelo	Calcium glycerophosph...	13007908	1
Caelo	Calcium glycerophosph...	13007912	2
Caelo	Calcium glycerophosph...	13A18-No2	1
Caelo	Calcium glycerophosph...	19100701	1
Caelo	Calcium glycerophosph...	15B03B03305050	1
Caelo	Calcium glycerophosph...	17359106	1
Caelo	Calcium glycerophosph...	13007902	1
Caelo	Calcium glycerophosph...	16231004	2
Caelo	Calcium glycerophosph...	14341009	1
Caelo	Calcium glycerophosph...	180918A	1
Caelo	Calcium glycerophosph...	17359107	3
Caelo	Calcium hydrogen phos...	12212816	1
Caelo	Calcium hydrogen phos...	14085301	1
Caelo	Calcium hydrogen phos...	12212813	1
Caelo	Calcium hydrogen phos...	14011401	1
Caelo	Calcium hydrogen phos...	10135228	1
Caelo	Calcium hydrogen phos...	12212804	5
Caelo	Calcium hydrogen phos...	12212807	1
Caelo	Calcium hydrogen phos...	10135216	1
Caelo	Calcium hydrogen phos...	11287501	1
Caelo	Calcium hydrogen phos...	12212806	1
Caelo	Calcium hydrogen phos...	13170701	5
Caelo	Calcium hydrogen phos...	12212822	1
Caelo	Calcium hydrogen phos...	11287509	2
Caelo	Calcium hydrogen phos...	12212811	5
Caelo	Calcium hydrogen phos...	12212812	2
Caelo	Calcium hydrogen phos...	10135201	1
Caelo	Calcium hydrogen phos...	13170704	2
Caelo	Calcium hydrogen phos...	11287510	1
Caelo	Calcium hydrogen phos...	101352009	1
Caelo	Calcium hydrogen phos...	10135223	1
Caelo	Calcium hydrogen phos...	11287503	2
Caelo	Calcium lactate penta...	16333111	4
Caelo	Calcium lactate penta...	17355503	1
Caelo	Calcium lactate penta...	17137707	2
Caelo	Calcium lactate penta...	16131312	1
Caelo	Calcium lactate penta...	16131302	1
Caelo	Calcium lactate penta...	18073413	2
Caelo	Calcium lactate penta...	15101412	3
Caelo	Calcium lactate penta...	16131307	1
Caelo	Calcium lactate penta...	15101409	1
Caelo	Calcium lactate penta...	16333106	3
Caelo	Calcium lactate penta...	13080802	2
Caelo	Calcium lactate penta...	14H01B05300502	1
Caelo	Calcium lactate penta...	11369606	1
Caelo	Calcium lactate penta...	12E10-N03	1
Caelo	Calcium lactate penta...	13080805	2
Caelo	Calcium lactate penta...	13080804	2
Caelo	Calcium lactate penta...	12C15-N02	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Calcium lactate penta...	11369610	1
Caelo	Calcium lactate penta...	15227804	2
Caelo	Calcium lactate penta...	11369618	1
Caelo	Calcium lactate penta...	11369614	2
Caelo	Calcium pantothenate	12340809	1
Caelo	Calcium pantothenate	12340832	1
Caelo	Calcium pantothenate	8041407	1
Caelo	Calcium pantothenate	12340846	1
Caelo	Calcium pantothenate	11105731	1
Caelo	Calcium pantothenate	12340808	3
Caelo	Calcium pantothenate	12236622	1
Caelo	Calcium pantothenate	12340859	1
Caelo	Calcium pantothenate	14208004	1
Caelo	Calcium pantothenate	15344207	1
Caelo	Calcium pantothenate	15344214	1
Caelo	Calcium pantothenate	17335208	3
Caelo	Calcium pantothenate	1013Q-03022	1
Caelo	Calcium pantothenate	27111303	1
Caelo	Calcium pantothenate	16251608	2
Caelo	Calcium pantothenate	12340819	2
Caelo	Calcium sulfate dihyd...	10132403	1
Caelo	Calcium sulfate dihyd...	12316301	1
Caelo	Calendula flower	12221510-3435311151	1
Caelo	Calendula flower	12221510-3356319071	1
Caelo	Calendula flower	122221707	1
Caelo	Calendula flower	12221710	1
Caelo	Calendula flower	12221705	1
Caelo	Calendula flower	12221511	1
Caelo	Calendula flower	12221508	1
Caelo	Calendula flower	13235502	3
Caelo	Calendula flower	15291501	1
Caelo	Calendula flower	12221708	3
Caelo	Calendula flower	13235503	3
Caelo	Calendula flower	12221706	2
Caelo	Calendula flower, cut	13410202	2
Caelo	Calendula flower, cut	12221705	4
Caelo	Calendula flower, cut	5021305	1
Caelo	Calendula flower, cut	12221701	1
Caelo	Calendula flower, cut	31121202	1
Caelo	Calendula flower, cut	11187611	3
Caelo	Calendula flower, cut	122221704	1
Caelo	Calendula flower, cut	1187612	2
Caelo	Calendula flower, cut	1325502	1
Caelo	Calendula flower, cut	27101201	1
Caelo	Calendula flower, cut	11187610	3
Caelo	Calendula flower, cut	11187612	11
Caelo	Calendula flower, cut	13235506	2
Caelo	Calendula flower, cut	13235505	1
Caelo	Calendula flower, cut	19111206	1
Caelo	Calendula flower, cut	21091301	1
Caelo	Calendula flower, cut	11221407	1
Caelo	Calendula flower, cut	18121301	1
Caelo	Calendula flower, cut	AKD1	1
Caelo	Calendula flower, cut	15040301	1
Caelo	Calendula flower, cut	12221708	11
Caelo	Calendula flower, cut	12221707	2
Caelo	Calendula flower, cut	295556	1
Caelo	Calendula flower, cut	12221710	2
Caelo	Calendula flower, cut	12221502	2
Caelo	Calendula flower, cut	13235502	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Calendula flower, cut	9112013C	1
Caelo	Calendula flower, cut	12071201	1
Caelo	Calendula flower, cut	13235501	4
Caelo	Calendula flower, cut	8041302	1
Caelo	Calendula flower, cut	12221704	11
Caelo	Calendula flower, cut	9121303	1
Caelo	Calendula flower, cut	3061301	1
Caelo	Calendula flower, cut	10291904	1
Caelo	Calendula flower, cut	22061203	1
Caelo	Camphor	13205517	1
Caelo	Camphor	13205522	1
Caelo	Camphor	18327905	1
Caelo	Camphor	13195316	6
Caelo	Camphor	12061411	1
Caelo	Camphor	15/18.B11-321658	1
Caelo	Camphor	13111708	1
Caelo	Camphor	13205522	2
Caelo	Camphor	13111707	1
Caelo	Camphor	16209310	1
Caelo	Camphor	15316823	2
Caelo	Camphor	13195314	5
Caelo	Camphor	3091203	1
Caelo	Camphor	13205502	1
Caelo	Camphor	2811Q-01221	1
Caelo	Camphor	11286905	2
Caelo	Camphor	13195318	1
Caelo	Camphor	17288307	1
Caelo	Camphor	13108118	2
Caelo	Camphor	4031503	1
Caelo	Camphor	18019806	9
Caelo	Camphor	14202801	3
Caelo	Camphor	16235404	2
Caelo	Camphor	14310109	1
Caelo	Camphor	14310715	1
Caelo	Camphor	13195322	7
Caelo	Camphor	6021901	1
Caelo	Camphor	914E-01221	1
Caelo	Camphor	14224303	1
Caelo	Camphor	62(17288306)	1
Caelo	Camphor	18108904	2
Caelo	Camphor	13K28-04-294937	1
Caelo	Camphor	14224303	3
Caelo	Camphor	18019811	2
Caelo	Camphor	11032023	1
Caelo	Camphor	15316805	1
Caelo	Camphor	17050606	2
Caelo	Camphor	162093029	1
Caelo	Camphor	13108106	2
Caelo	Camphor	4101303	1
Caelo	Camphor	13195314	3
Caelo	Camphor	12054151	2
Caelo	Camphor	12054149	9
Caelo	Camphor	12054114	2
Caelo	Camphor	13205501	2
Caelo	Camphor	13195304	1
Caelo	Camphor	12054143	1
Caelo	Camphor	13195301	2
Caelo	Camphor	13195311	1
Caelo	Camphor	12243701	1
Caelo	Camphor	713A-01221	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Camphor	13205512	1
Caelo	Camphor	13205511	1
Caelo	Camphor	11286910	1
Caelo	Camphor	13205511	2
Caelo	Camphor	2344905	1
Caelo	Camphor	12228205	1
Caelo	Camphor	11074509	1
Caelo	Camphor	17004304	1
Caelo	Camphor	11286909	2
Caelo	Camphor	12054105	1
Caelo	Camphor	12054107	4
Caelo	Camphor	13195311	1
Caelo	Camphor	12054124	2
Caelo	Camphor	25041401	1
Caelo	Camphor	12054113	1
Caelo	Camphor	120564113	1
Caelo	Camphor	17004302	3
Caelo	Camphor	12054141	3
Caelo	Camphor	12344905	2
Caelo	Camphor	17050605	14
Caelo	Camphor	12054142	6
Caelo	Camphor	12228206	1
Caelo	Camphor	18198003	1
Caelo	Camphor	17004301	14
Caelo	Camphor		1
Caelo	Camphor	17324104	4
Caelo	Camphor	14310710	2
Caelo	Camphor	16209304	3
Caelo	Camphor	18019803	4
Caelo	Camphor	15086806	1
Caelo	Camphor	15120902	1
Caelo	Camphor	1412023-01	5
Caelo	Camphor	15062702	3
Caelo	Camphor	15097802	2
Caelo	Camphor	15086801	1
Caelo	Camphor	14310705	1
Caelo	Camphor	15316801	4
Caelo	Camphor	14310705	13
Caelo	Camphor	315E-01221	1
Caelo	Camphor	5315I-01221	1
Caelo	Camphor	16209301	8
Caelo	Camphor	17288301	5
Caelo	Camphor	1410710	1
Caelo	Camphor	15316806	2
Caelo	Camphor	15320604	4
Caelo	Camphor	14310713	9
Caelo	Camphor	15397705	1
Caelo	Camphor	4444A150721	1
Caelo	Camphor	15120901	2
Caelo	Camphor	14310709	5
Caelo	Camphor	15316815	4
Caelo	Camphor	15316812	1
Caelo	Camphor	150622C	1
Caelo	Camphor	15316807	6
Caelo	Camphor	140918C	1
Caelo	Camphor	17134703	1
Caelo	Camphor	15316819	6
Caelo	Camphor	13195323	1
Caelo	Camphor	15397704	5
Caelo	Camphor	15397701	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Camphor	7071801	1
Caelo	Camphor	1531607	1
Caelo	Camphor	444A170221	1
Caelo	Camphor	16209309	8
Caelo	Camphor	14310722	2
Caelo	Camphor	17050605/1140418	1
Caelo	Camphor	17288305	1
Caelo	Camphor	14202807	3
Caelo	Camphor	16209306	1
Caelo	Camphor	14310703	1
Caelo	Camphor	7061802	1
Caelo	Camphor	16235405	2
Caelo	Camphor	16209305	6
Caelo	Camphor	15097803	2
Caelo	Camphor	13195326	3
Caelo	Camphor	14224304	1
Caelo	Camphor	14310718	4
Caelo	Camphor	13205512	1
Caelo	Camphor	15384201	1
Caelo	Camphor	12091805	1
Caelo	Camphor	2803020-01	1
Caelo	Camphor	12091808	1
Caelo	Cannabidiol	CBD-API1803	1
Caelo	Cannabidiol	110216-02	3
Caelo	Cannabidiol	CBD-API17-006	6
Caelo	Cannabidiol	CBD-API1801	5
Caelo	Captopril	14196905	3
Caelo	Captopril	15189510	2
Caelo	Captopril	14196901	1
Caelo	Captopril	15189514	6
Caelo	Captopril	14196904	8
Caelo	Captopril	5102-14-042	1
Caelo	Captopril	15189509	2
Caelo	Captopril	15189511	4
Caelo	Captopril	14B26-B01-310127	1
Caelo	Captopril	12F29-N04	1
Caelo	Captopril	12339605	4
Caelo	Captopril	15189502	1
Caelo	Captopril	15189501	6
Caelo	Captopril	15189517	5
Caelo	Captopril	14B26-B40-323267	4
Caelo	Captopril	12339604	2
Caelo	Captopril	12339610	2
Caelo	Captopril	14196903	1
Caelo	Captopril	12339602	1
Caelo	Captopril	15189503	1
Caelo	Captopril	15189515	1
Caelo	Captopril	12339619	1
Caelo	Captopril	12339616	2
Caelo	Captopril	14196903	1
Caelo	Captopril	15189512	6
Caelo	Captopril	19011802	2
Caelo	Captopril	17077208	8
Caelo	Captopril	15189526	2
Caelo	Captopril	17210205	9
Caelo	Captopril	18058009	1
Caelo	Captopril	15189524	2
Caelo	Captopril	1519521	1
Caelo	Captopril	17077205	1
Caelo	Captopril	15189525	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Captopril	17077202	6
Caelo	Captopril	17077209	13
Caelo	Captopril	1805005	1
Caelo	Captopril	18058005	1
Caelo	Captopril	17077206	1
Caelo	Captopril	17077201	3
Caelo	Captopril	18058006	2
Caelo	Captopril	18058002	2
Caelo	Captopril	18058003	1
Caelo	Captopril	18058004	2
Caelo	Captopril	17210202	8
Caelo	Captopril	15189521	5
Caelo	Captopril	17G27-B02	2
Caelo	Captopril	15189519	2
Caelo	Captopril	16B26-B01310127	1
Caelo	Captopril	19011801	1
Caelo	Caraway fruit	13091401	1
Caelo	Caraway fruit	12226604	5
Caelo	Caraway fruit	13397104	2
Caelo	Caraway fruit	12226603	7
Caelo	Caraway fruit	12226601	1
Caelo	Caraway fruit	12226611	3
Caelo	Caraway fruit	12226608	1
Caelo	Caraway fruit	13306004	1
Caelo	Caraway fruit	13397105	8
Caelo	Caraway fruit	811E-01537	1
Caelo	Caraway fruit	15097101	8
Caelo	Caraway fruit	1509101	1
Caelo	Caraway fruit	15097116	1
Caelo	Caraway fruit	13306008	1
Caelo	Caraway fruit	15409103	2
Caelo	Caraway fruit	15097108	4
Caelo	Caraway fruit	16308502	1
Caelo	Caraway fruit	12226612	8
Caelo	Caraway fruit	15097107	2
Caelo	Caraway fruit	15097111	1
Caelo	Caraway fruit	13041903	1
Caelo	Caraway fruit	15097115	6
Caelo	Caraway fruit	12343403	1
Caelo	Caraway fruit	13306006	4
Caelo	Caraway fruit	13306005	2
Caelo	Caraway fruit	13306002	5
Caelo	Caraway fruit	11003211	1
Caelo	Caraway fruit	2103311537	1
Caelo	Caraway fruit	9267815	1
Caelo	Caraway fruit	13397102	2
Caelo	Caraway fruit	11003207	1
Caelo	Caraway fruit	413M	1
Caelo	Caraway fruit	15409105	1
Caelo	Caraway fruit	17256303	1
Caelo	Caraway fruit	16308501	1
Caelo	Caraway fruit	16308503	1
Caelo	Caraway fruit	18121717	1
Caelo	Caraway fruit	12226608	4
Caelo	Caraway fruit	12226504	1
Caelo	Caraway fruit	12226608	1
Caelo	Caraway fruit	12226606	2
Caelo	Caraway fruit	2041017/16308501	1
Caelo	Caraway fruit	16165308	1
Caelo	Caraway fruit	290416	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Caraway fruit	12226609	5
Caelo	Caraway fruit	11121704	1
Caelo	Caraway fruit	4061303	1
Caelo	Caraway fruit	11003217	5
Caelo	Caraway fruit	15097103	3
Caelo	Caraway fruit	13041901	1
Caelo	Caraway fruit	12226614	2
Caelo	Caraway fruit	11378309	1
Caelo	Caraway fruit	12226501	5
Caelo	Caraway fruit	15097109	1
Caelo	Caraway fruit	27121202	1
Caelo	Caraway fruit	13306003	1
Caelo	Caraway fruit	15409104	2
Caelo	Caraway fruit, pounded	12116904	3
Caelo	Caraway fruit, pounded	13367101	1
Caelo	Carbomers 35000 / 50000	1510801	1
Caelo	Carbomers 35000 / 50000	15101810	5
Caelo	Carbomers 35000 / 50000	15101807	3
Caelo	Carbomers 35000 / 50000	181016184	1
Caelo	Carbomers 35000 / 50000	162375	1
Caelo	Carbomers 35000 / 50000	18307403	2
Caelo	Carbomers 35000 / 50000	18092412	1
Caelo	Carbomers 35000 / 50000	18307406	1
Caelo	Carbomers 35000 / 50000	18307401	1
Caelo	Carbomers 35000 / 50000	15101801	5
Caelo	Carbomers 35000 / 50000	15101802	5
Caelo	Carbomers 35000 / 50000	13186207	1
Caelo	Carbomers 35000 / 50000	15015703	1
Caelo	Carbomers 35000 / 50000	15015702	1
Caelo	Carbomers 35000 / 50000	15101809	4
Caelo	Carbomers 35000 / 50000	31631293	1
Caelo	Carbomers 35000 / 50000	15015709	2
Caelo	Carbomers 35000 / 50000	15439002	6
Caelo	Carbomers 35000 / 50000	162530	1
Caelo	Carbomers 35000 / 50000	16253014	3
Caelo	Carbomers 35000 / 50000	18092401	12
Caelo	Carbomers 35000 / 50000	17146501	2
Caelo	Carbomers 35000 / 50000	17361403	2
Caelo	Carbomers 35000 / 50000	17053182	1
Caelo	Carbomers 35000 / 50000	18092402	4
Caelo	Carbomers 35000 / 50000	16253012	4
Caelo	Carbomers 35000 / 50000	1081702	1
Caelo	Carbomers 35000 / 50000	18049304	1
Caelo	Carbomers 35000 / 50000	11676321	1
Caelo	Carbomers 35000 / 50000	16253002	2
Caelo	Carbomers 35000 / 50000	16253009	8
Caelo	Carbomers 35000 / 50000	18092411	1
Caelo	Carbomers 35000 / 50000	16237501	2
Caelo	Carbomers 35000 / 50000	16253001	1
Caelo	Carbomers 35000 / 50000	16253005	4
Caelo	Carbomers 35000 / 50000	16253004	9
Caelo	Carbomers 35000 / 50000	18092405	5
Caelo	Carbomers 35000 / 50000	18092406	1
Caelo	Carbomers 35000 / 50000	9011903	1
Caelo	Carbomers 35000 / 50000	154390002	1
Caelo	Carbomers 35000 / 50000	15439003	11
Caelo	Carbomers 35000 / 50000	15257401	1
Caelo	Carbomers 35000 / 50000	15257403	1
Caelo	Carbomers 35000 / 50000	2215M-01232	1
Caelo	Carbomers 35000 / 50000	17361402	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Carbomers 35000 / 50000	18092408	2
Caelo	Carbomers 35000 / 50000	16253007	4
Caelo	Carbomers 35000 / 50000	14149910	2
Caelo	Carbomers 35000 / 50000	2913M1232	1
Caelo	Carbomers 35000 / 50000	13194507	3
Caelo	Carbomers 35000 / 50000	14057902	2
Caelo	Carbomers 35000 / 50000	12274908	1
Caelo	Carbomers 35000 / 50000	14149911	2
Caelo	Carbomers 35000 / 50000	12274908	3
Caelo	Carbomers 35000 / 50000	13194503	1
Caelo	Carbomers 35000 / 50000	13194512	1
Caelo	Carbomers 35000 / 50000	11307511	1
Caelo	Carbomers 35000 / 50000	14016208	7
Caelo	Carbomers 35000 / 50000	13194511	1
Caelo	Carbomers 35000 / 50000	14057903	2
Caelo	Carbomers 35000 / 50000	12274907	6
Caelo	Carbomers 35000 / 50000	11351605	2
Caelo	Carbomers 35000 / 50000	11351623	1
Caelo	Carbomers 35000 / 50000	11351618	2
Caelo	Carbomers 35000 / 50000	12274907	1
Caelo	Carbomers 35000 / 50000	11013605	1
Caelo	Carbomers 35000 / 50000	13186201	2
Caelo	Carbomers 35000 / 50000	13194501	4
Caelo	Carbomers 35000 / 50000	13194501	1
Caelo	Carbomers 35000 / 50000	13194503	8
Caelo	Carbomers 35000 / 50000	11351604	1
Caelo	Carbomers 35000 / 50000	12274912	4
Caelo	Carbomers 35000 / 50000	13186203	3
Caelo	Carbomers 35000 / 50000	14149906	3
Caelo	Carbomers 35000 / 50000	28111310	1
Caelo	Carbomers 35000 / 50000	12274914	3
Caelo	Carbomers 35000 / 50000	12274916	9
Caelo	Carbomers 35000 / 50000	13194507	1
Caelo	Carbomers 35000 / 50000	14016209	1
Caelo	Carbomers 35000 / 50000	11013605	1
Caelo	Carbomers 35000 / 50000	13186205	2
Caelo	Carbomers 35000 / 50000	12274910	1
Caelo	Carbomers 35000 / 50000	11351618	1
Caelo	Carbomers 35000 / 50000	14149904	1
Caelo	Carbomers 35000 / 50000	13186205	4
Caelo	Cardamom fruit	13026003	1
Caelo	Cardamom fruit	13260502-3537411072	1
Caelo	Cardamom fruit	11309110	1
Caelo	Carmellose sodium (40...	14149806	1
Caelo	Carmellose sodium (40...	15123908	3
Caelo	Carmellose sodium (40...	12354416	2
Caelo	Carmellose sodium (40...	13104813	2
Caelo	Carmellose sodium (40...	28051402	1
Caelo	Carmellose sodium (40...	13299302	1
Caelo	Carmellose sodium (40...	17348605	1
Caelo	Carmellose sodium (40...	13G19-N02	1
Caelo	Carmellose sodium (40...	15123901	3
Caelo	Carmellose sodium (40...	14029211	1
Caelo	Carmellose sodium (40...	13104802	2
Caelo	Carmellose sodium (40...	13104807	1
Caelo	Carmellose sodium (40...	13104802	1
Caelo	Carmellose sodium (40...	16171402	3
Caelo	Carmellose sodium (40...	15123909	1
Caelo	Carmellose sodium (40...	12354408	1
Caelo	Carmellose sodium (40...	11124211	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Carmellose sodium (40...	12255708	1
Caelo	Carmellose sodium (40...	13104810	3
Caelo	Carmellose sodium (40...	25041805	1
Caelo	Carmellose sodium (40...	17062309	2
Caelo	Carmellose sodium (40...	15123912	1
Caelo	Carmellose sodium (40...	12354414	1
Caelo	Carmellose sodium (40...	17348601	2
Caelo	Carmellose sodium (40...	17062303	1
Caelo	Carmellose sodium (40...	17052304	1
Caelo	Carmellose sodium (40...	13299302	1
Caelo	Carmellose sodium (40...	17062304	4
Caelo	Carmellose sodium (40...	17348606	1
Caelo	Carmellose sodium (40...	13104807	1
Caelo	Carmellose sodium (40...	14029201	1
Caelo	Carmellose sodium (40...	1315Q-03302	1
Caelo	Carmellose sodium (40...	12354416	1
Caelo	Carmellose sodium (40...	17062307	4
Caelo	Cayenne pepper, powde...	13205106	1
Caelo	Celery root, cut	12319108	1
Caelo	Celery root, cut	12319107	1
Caelo	Celery root, cut	10270012	1
Caelo	Celery root, cut	12319105	2
Caelo	Celery root, cut	10270009	1
Caelo	Cellulose, microcryst...	160402901	1
Caelo	Centaury, cut	12268601	3
Caelo	Centaury, cut	12268602	19
Caelo	Centaury, cut	8021302	1
Caelo	Centaury, cut	10104012	1
Caelo	Centaury, cut	11097609	3
Caelo	Centaury, cut	11097605	2
Caelo	Centaury, cut	3711B120628	1
Caelo	Centaury, cut	12268605	4
Caelo	Centaury, cut	15224102	1
Caelo	Centaury, cut	11097602	1
Caelo	Centaury, cut	12268608	10
Caelo	Centaury, cut	12379601	1
Caelo	Centaury, cut	12268606	13
Caelo	Cetomacrogol emulsify...	16348101	1
Caelo	Cetomacrogol emulsify...	16080403	1
Caelo	Cetomacrogol emulsify...	16080407	1
Caelo	Cetomacrogol emulsify...	15277505	1
Caelo	Cetomacrogol emulsify...	15277501	1
Caelo	Cetomacrogol emulsify...	11001603Nachprüfung	1
Caelo	Cetomacrogol emulsify...	15113703	1
Caelo	Cetomacrogol emulsify...	14055002	1
Caelo	Cetostearyl alcohol	13206206	4
Caelo	Cetostearyl alcohol	14054504	1
Caelo	Cetostearyl alcohol	14054508	2
Caelo	Cetostearyl alcohol	15272709	5
Caelo	Cetostearyl alcohol	15272702	1
Caelo	Cetostearyl alcohol	15271405	2
Caelo	Cetostearyl alcohol	15090701	2
Caelo	Cetostearyl alcohol	14158201	1
Caelo	Cetostearyl alcohol	14215007	6
Caelo	Cetostearyl alcohol	14346201	4
Caelo	Cetostearyl alcohol	14322506	1
Caelo	Cetostearyl alcohol	18003404	3
Caelo	Cetostearyl alcohol	13018306	1
Caelo	Cetostearyl alcohol	13018305	4
Caelo	Cetostearyl alcohol	17049209	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cetostearyl alcohol	14061317	1
Caelo	Cetostearyl alcohol	12065709	2
Caelo	Cetostearyl alcohol	12021304	1
Caelo	Cetostearyl alcohol	13233002	1
Caelo	Cetostearyl alcohol	14158201	1
Caelo	Cetostearyl alcohol	14054508	1
Caelo	Cetostearyl alcohol	14054504	3
Caelo	Cetostearyl alcohol	18018001	1
Caelo	Cetostearyl alcohol	17321603	1
Caelo	Cetostearyl alcohol	18A23-B04	1
Caelo	Cetostearyl alcohol	17139201	1
Caelo	Cetostearyl alcohol	16238607	1
Caelo	Cetostearyl alcohol	16238604	1
Caelo	Cetostearyl alcohol	160914	1
Caelo	Cetostearyl alcohol	171215CSA/15272717	1
Caelo	Cetostearyl alcohol	13-B01-327355	1
Caelo	Cetostearyl alcohol	17198003	2
Caelo	Cetostearyl alcohol	8041708	1
Caelo	Cetostearyl alcohol	15272717	5
Caelo	Cetostearyl alcohol	16119604	3
Caelo	Cetostearyl alcohol	15272709	1
Caelo	Cetostearyl alcohol	17321303	3
Caelo	Cetostearyl alcohol	1314405	1
Caelo	Cetostearyl alcohol	13144105	1
Caelo	Cetostearyl alcohol	14346202	1
Caelo	Cetostearyl alcohol	17132001	1
Caelo	Cetostearyl alcohol	16048812	1
Caelo	Cetostearyl alcohol	15271407	1
Caelo	Cetostearyl alcohol	18018005	2
Caelo	Cetostearyl alcohol	17261702	1
Caelo	Cetostearyl alcohol	11162305	1
Caelo	Cetostearyl alcohol	12025301	2
Caelo	Cetyl alcohol	16335802	7
Caelo	Cetyl alcohol	16335805	1
Caelo	Cetyl alcohol	11153405	1
Caelo	Cetyl alcohol	17047307	6
Caelo	Cetyl alcohol	12271103	5
Caelo	Cetyl alcohol	28101701	1
Caelo	Cetyl alcohol	17047306	3
Caelo	Cetyl alcohol	17234902	1
Caelo	Cetyl alcohol	12338604	6
Caelo	Cetyl alcohol	1212E-01076	1
Caelo	Cetyl alcohol	13096901	6
Caelo	Cetyl alcohol	20250112-1	1
Caelo	Cetyl alcohol	17234903	1
Caelo	Cetyl alcohol	12018605	2
Caelo	Cetyl alcohol	11320702	2
Caelo	Cetyl alcohol	7071807	1
Caelo	Cetyl alcohol	17321612	6
Caelo	Cetyl alcohol	18153102	7
Caelo	Cetyl alcohol	15323316	3
Caelo	Cetyl alcohol	1112013B	1
Caelo	Cetyl alcohol	13096904	2
Caelo	Cetyl alcohol	17321603	6
Caelo	Cetyl alcohol	14262204	3
Caelo	Cetyl alcohol	16335807	6
Caelo	Cetyl alcohol	13286502	4
Caelo	Cetyl alcohol	15323311	7
Caelo	Cetyl alcohol	15323301	1
Caelo	Cetyl alcohol	13210504	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cetyl alcohol	13210504	1
Caelo	Cetyl alcohol	13286502	1
Caelo	Cetyl alcohol	16106307	3
Caelo	Cetyl alcohol	15323301	6
Caelo	Cetyl alcohol	14262204	2
Caelo	Cetyl alcohol	15323308	1
Caelo	Cetyl alcohol	14349101	14
Caelo	Cetyl alcohol	1523316	1
Caelo	Cetyl alcohol	15323311	1
Caelo	Cetyl palmitate 15	17164004	9
Caelo	Cetyl palmitate 15	14351606	4
Caelo	Cetyl palmitate 15	17043504	11
Caelo	Cetyl palmitate 15	14213706	1
Caelo	Cetyl palmitate 15	14351612	15
Caelo	Cetyl palmitate 15	ch.b.15301705	1
Caelo	Cetyl palmitate 15	17164005	17
Caelo	Cetyl palmitate 15	14351609	2
Caelo	Cetyl palmitate 15	14351602	5
Caelo	Cetyl palmitate 15	14041505	1
Caelo	Cetyl palmitate 15	14213710	1
Caelo	Cetyl palmitate 15	14213709	1
Caelo	Cetyl palmitate 15	17271506	1
Caelo	Cetyl palmitate 15	11021502	1
Caelo	Cetyl palmitate 15	14060709	2
Caelo	Cetyl palmitate 15	14213709	3
Caelo	Cetyl palmitate 15	17164007	1
Caelo	Cetyl palmitate 15	17271505	6
Caelo	Cetyl palmitate 15	17271504	5
Caelo	Cetyl palmitate 15	14351601	2
Caelo	Cetyl palmitate 15	14351612	1
Caelo	Cetyl palmitate 15	142137710	1
Caelo	Cetyl palmitate 15	17043503	4
Caelo	Cetyl palmitate 15	16053804	1
Caelo	Cetyl palmitate 15	13131905Nachprüfen	1
Caelo	Cetyl palmitate 15	16053805	2
Caelo	Cetyl palmitate 15	15301711	7
Caelo	Cetyl palmitate 15	15301708	1
Caelo	Cetyl palmitate 15	16053803	9
Caelo	Cetyl palmitate 15	15301701	1
Caelo	Cetyl palmitate 15	15301705	8
Caelo	Cetyl palmitate 15	15632	1
Caelo	Cetyl palmitate 15	14351607	1
Caelo	Cetyl palmitate 15	15301709	6
Caelo	Cetyl palmitate 15	15301704	5
Caelo	Cetyl palmitate 15	17271502	5
Caelo	Cetyl palmitate 15	14351614	1
Caelo	Cetyl palmitate 15	15301701	7
Caelo	Cetyl palmitate 15	16053811	8
Caelo	Cetyl palmitate 15	14351609	1
Caelo	Cetyl palmitate 15	16186805	20
Caelo	Cetyl palmitate 15	16186809	7
Caelo	Cetyl palmitate 15	28041731	1
Caelo	Cetyl palmitate 15	16186806	1
Caelo	Cetyl palmitate 15	15301703	8
Caelo	Cetyl palmitate 15	14060705	1
Caelo	Cetyl palmitate 15	17043509	1
Caelo	Cetyl palmitate 15	12018204/Int212191	1
Caelo	Cetyl palmitate 15	26071201	1
Caelo	Cetyl palmitate 15	12082901	1
Caelo	Cetyl palmitate 15	12018208	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cetyl palmitate 15	12018215	1
Caelo	Cetyl palmitate 15	14060705	6
Caelo	Cetyl palmitate 15	11042019	1
Caelo	Cetyl palmitate 15	12228808	4
Caelo	Cetyl palmitate 15	12018204	3
Caelo	Cetyl palmitate 15	12228805	1
Caelo	Cetyl palmitate 15	19121803	1
Caelo	Cetyl palmitate 15	1228808	1
Caelo	Cetyl palmitate 15	18273704	2
Caelo	Cetyl palmitate 15	13295910	8
Caelo	Cetyl palmitate 15	1313190	1
Caelo	Cetyl palmitate 15	30101301	1
Caelo	Cetyl palmitate 15	10068410	1
Caelo	Cetyl palmitate 15	13295919	11
Caelo	Cetyl palmitate 15	14060709	4
Caelo	Cetyl palmitate 15	13131914	4
Caelo	Cetyl palmitate 15	18149902	1
Caelo	Cetyl palmitate 15	13295921	1
Caelo	Cetyl palmitate 15	13131911	2
Caelo	Cetyl palmitate 15	14060706	4
Caelo	Cetyl palmitate 15	12331916	8
Caelo	Cetyl palmitate 15	12331919	2
Caelo	Cetyl palmitate 15	13131906	8
Caelo	Cetyl palmitate 15	12228818	5
Caelo	Cetyl palmitate 15	13131905	4
Caelo	Cetyl palmitate 15	14213705	3
Caelo	Cetyl palmitate 15	14060706	2
Caelo	Cetyl palmitate 15	12228801/Int212191	1
Caelo	Ceylonese Cinnamon, p...	12242818	1
Caelo	Ceylonese Cinnamon, p...	12242806	2
Caelo	Ceylonese Cinnamon, p...	12242815	1
Caelo	Chelidonium majus her...	311092	1
Caelo	Chickweed, cut	13040501	1
Caelo	Chickweed, cut	4551B130902	1
Caelo	Chickweed, cut	13254404	1
Caelo	Chickweed, cut	13040502	2
Caelo	Chickweed, cut	13254401	1
Caelo	Chickweed, cut	13040503	1
Caelo	Chickweed, cut	13040504	1
Caelo	Chickweed, cut	13254403	2
Caelo	Chickweed, cut	13040502-3283310662	1
Caelo	Chickweed, cut	14212411	1
Caelo	Chicory herb, cut	12369102	3
Caelo	Chicory herb, cut	12354306	1
Caelo	Chicory herb, cut	12354302	3
Caelo	Chicory herb, cut	14031315	1
Caelo	Chicory herb, cut	12354304	1
Caelo	Chicory herb, cut	12034302	1
Caelo	Chicory herb, cut	12034301	2
Caelo	Chicory herb, cut	12034303	3
Caelo	Chicory root, cut	11374504	3
Caelo	Chicory root, cut	6221As101202	1
Caelo	Chicory root, cut	11081803	1
Caelo	Chicory root, cut	11374510	2
Caelo	Chicory root, cut	6221A101202	4
Caelo	Chicory root, cut	11374506	1
Caelo	Chicory root, cut	11374501	1
Caelo	Chicory root, cut	91338319	1
Caelo	Chloral hydrate	16248703	1
Caelo	Chloral hydrate	15448313	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chloral hydrate	17346302	10
Caelo	Chloral hydrate	15448307	5
Caelo	Chloral hydrate	15448315	5
Caelo	Chloral hydrate	17346309	3
Caelo	Chloral hydrate	12017518	1
Caelo	Chloral hydrate	16355811	1
Caelo	Chloral hydrate	12017503	1
Caelo	Chloral hydrate	12017512	1
Caelo	Chloral hydrate	13002425	4
Caelo	Chloral hydrate	17346303	4
Caelo	Chloral hydrate	13d22-n11	1
Caelo	Chloral hydrate	25021401	1
Caelo	Chloral hydrate	15448314	7
Caelo	Chloral hydrate	13213311	1
Caelo	Chloral hydrate	17346311	1
Caelo	Chloral hydrate	14038301	2
Caelo	Chloral hydrate	17346306	18
Caelo	Chloral hydrate	14038311	4
Caelo	Chloral hydrate	173463021	1
Caelo	Chloral hydrate	17346307	4
Caelo	Chloral hydrate	17314312	1
Caelo	Chloral hydrate	30081807	1
Caelo	Chloral hydrate	15448312	10
Caelo	Chloral hydrate	14101307	1
Caelo	Chloral hydrate	15166507	10
Caelo	Chloral hydrate	17314322	1
Caelo	Chloral hydrate	15448304	2
Caelo	Chloral hydrate	17314318	1
Caelo	Chloral hydrate	17121402	1
Caelo	Chloral hydrate	6051703	1
Caelo	Chloral hydrate	15448318	7
Caelo	Chloral hydrate	14218714	4
Caelo	Chloral hydrate	15448310	1
Caelo	Chloral hydrate	14218712	2
Caelo	Chloral hydrate	14218720	1
Caelo	Chloral hydrate	160906C	1
Caelo	Chloral hydrate	14218719	4
Caelo	Chloral hydrate	14218713	1
Caelo	Chloral hydrate	17550601	1
Caelo	Chloral hydrate	15083703	8
Caelo	Chloral hydrate	15166503	4
Caelo	Chloral hydrate	15083705	1
Caelo	Chloral hydrate	15166501	3
Caelo	Chloral hydrate	15166505	6
Caelo	Chloral hydrate	15448303	3
Caelo	Chloral hydrate	15166502	3
Caelo	Chloral hydrate	15166506	6
Caelo	Chloral hydrate	15163901	1
Caelo	Chloral hydrate	15166508	5
Caelo	Chloral hydrate	15166513	2
Caelo	Chloral hydrate	15448306	7
Caelo	Chloral hydrate	15083704	1
Caelo	Chloral hydrate	14218709	5
Caelo	Chloral hydrate	3121403	1
Caelo	Chloral hydrate	14218705	4
Caelo	Chloral hydrate	15448325	2
Caelo	Chloral hydrate	17346310	2
Caelo	Chloral hydrate	16355804	2
Caelo	Chloral hydrate	19091408	1
Caelo	Chloral hydrate	17247702	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chloral hydrate	15448320	4
Caelo	Chloral hydrate	16355803	1
Caelo	Chloral hydrate	15448324	6
Caelo	Chloral hydrate	15448326	7
Caelo	Chloral hydrate	14038306	1
Caelo	Chloral hydrate	17314301	3
Caelo	Chloral hydrate	1403815	1
Caelo	Chloral hydrate	17520601	1
Caelo	Chloral hydrate	14218704	3
Caelo	Chloral hydrate	15448322	5
Caelo	Chloral hydrate	14218701	1
Caelo	Chloral hydrate	17247704	5
Caelo	Chloral hydrate	15448321	4
Caelo	Chloramphenicol	4712Q-01259	1
Caelo	Chloramphenicol	1137070	1
Caelo	Chloramphenicol	11370706	3
Caelo	Chloramphenicol	12G02-N01	12
Caelo	Chloramphenicol	12072606	1
Caelo	Chloramphenicol	10341703	1
Caelo	Chloramphenicol	11370703	2
Caelo	Chloramphenicol	11370705	3
Caelo	Chloramphenicol	11370707	2
Caelo	Chloramphenicol	10341717	1
Caelo	Chloramphenicol	11370702	3
Caelo	Chloramphenicol	14246917	1
Caelo	Chloramphenicol	10341711	2
Caelo	Chloramphenicol	15358903	3
Caelo	Chloramphenicol	15358905	19
Caelo	Chloramphenicol	2181116	1
Caelo	Chloramphenicol	15358907	3
Caelo	Chloramphenicol	150528C	1
Caelo	Chloramphenicol	14346914	3
Caelo	Chloramphenicol	14346922	1
Caelo	Chloramphenicol	3081602	1
Caelo	Chloramphenicol	14346918	17
Caelo	Chloramphenicol	1535896	1
Caelo	Chloramphenicol	15358909	6
Caelo	Chloramphenicol	130416-2	1
Caelo	Chloramphenicol	14356910	1
Caelo	Chloramphenicol	14346917	7
Caelo	Chloramphenicol	14346911	20
Caelo	Chloramphenicol	14346913	8
Caelo	Chloramphenicol	1436908	1
Caelo	Chloramphenicol	1434608	1
Caelo	Chloramphenicol	14346912	4
Caelo	Chloramphenicol	14346908	9
Caelo	Chloramphenicol	14346903	1
Caelo	Chloramphenicol	14346907	8
Caelo	Chloramphenicol	130416-1	1
Caelo	Chloramphenicol	13080925	11
Caelo	Chloramphenicol	15358906	11
Caelo	Chloramphenicol	16278609	12
Caelo	Chloramphenicol	18017814	1
Caelo	Chloramphenicol	8017811	1
Caelo	Chloramphenicol	18017811	4
Caelo	Chloramphenicol	181017809	1
Caelo	Chloramphenicol	18017809	10
Caelo	Chloramphenicol	18017804	2
Caelo	Chloramphenicol	18017806	11
Caelo	Chloramphenicol	1704021-03	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chloramphenicol	18017802	2
Caelo	Chloramphenicol	17343803	14
Caelo	Chloramphenicol	15358908	10
Caelo	Chloramphenicol	16278611	1
Caelo	Chloramphenicol	17343802	4
Caelo	Chloramphenicol	16258612	1
Caelo	Chloramphenicol	16358909	1
Caelo	Chloramphenicol	15358912	7
Caelo	Chloramphenicol	16278609/180420CA	1
Caelo	Chloramphenicol	160200201	1
Caelo	Chloramphenicol	16278612	12
Caelo	Chloramphenicol	16278009	1
Caelo	Chloramphenicol	16278610	1
Caelo	Chloramphenicol		1
Caelo	Chloramphenicol	1627862	1
Caelo	Chloramphenicol	14346904	4
Caelo	Chloramphenicol	130809	1
Caelo	Chloramphenicol	14346901	4
Caelo	Chloramphenicol	13A31-B02	9
Caelo	Chloramphenicol	16278608	9
Caelo	Chloramphenicol	310311	1
Caelo	Chloramphenicol	11370713	3
Caelo	Chloramphenicol	16278602	3
Caelo	Chloramphenicol	1627898	1
Caelo	Chloramphenicol	16278606	6
Caelo	Chloramphenicol	12248509	14
Caelo	Chloramphenicol	DY00911233	4
Caelo	Chloramphenicol	12248505	2
Caelo	Chloramphenicol	1137004	1
Caelo	Chloramphenicol	12248507	10
Caelo	Chloramphenicol	12248506	8
Caelo	Chloramphenicol	14346906	3
Caelo	Chloramphenicol	1137714	1
Caelo	Chloramphenicol	11370715	1
Caelo	Chloramphenicol	11370714	2
Caelo	Chloramphenicol	11370717	2
Caelo	Chloramphenicol	3125204	1
Caelo	Chloramphenicol	12248501	15
Caelo	Chloramphenicol	12g02n01	1
Caelo	Chloramphenicol	10341707	1
Caelo	Chloramphenicol	17A23-B08-337138	1
Caelo	Chloramphenicol	15358914	16
Caelo	Chloramphenicol	11370716	4
Caelo	Chloramphenicol	13080902	1
Caelo	Chloramphenicol	93323801	1
Caelo	Chloramphenicol	16/12-Bo7-329313	1
Caelo	Chloramphenicol	153589	1
Caelo	Chloramphenicol	13080924	15
Caelo	Chloramphenicol	4412Q-01259	1
Caelo	Chloramphenicol	13080918	2
Caelo	Chloramphenicol	13080917	4
Caelo	Chloramphenicol	13080923	11
Caelo	Chloramphenicol	16278601	1
Caelo	Chloramphenicol	140704301	6
Caelo	Chloramphenicol	DY130525	2
Caelo	Chloramphenicol	13080907	14
Caelo	Chloramphenicol	13080915	1
Caelo	Chloramphenicol	15358913	22
Caelo	Chloramphenicol	12248525	6
Caelo	Chloramphenicol	13080903	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chloramphenicol	13080904	6
Caelo	Chloramphenicol	28011401	1
Caelo	Chloramphenicol	13080303	1
Caelo	Chloramphenicol	20102012J1	1
Caelo	Chloramphenicol	12248520	7
Caelo	Chloramphenicol	12248524	7
Caelo	Chloramphenicol	11031402	1
Caelo	Chlorhexidine diacetate	15261012	8
Caelo	Chlorhexidine diacetate	15261009	12
Caelo	Chlorhexidine diacetate	31032014B	1
Caelo	Chlorhexidine diacetate	15261004	2
Caelo	Chlorhexidine diacetate	15261003	7
Caelo	Chlorhexidine diacetate	15261007	1
Caelo	Chlorhexidine diacetate	15261011	1
Caelo	Chlorhexidine diacetate	15438903	1
Caelo	Chlorhexidine diacetate	15261014	22
Caelo	Chlorhexidine diacetate	15261007	13
Caelo	Chlorhexidine diacetate	15261009	1
Caelo	Chlorhexidine diacetate	12142734	1
Caelo	Chlorhexidine diacetate	13177009	1
Caelo	Chlorhexidine diacetate	23012015D	1
Caelo	Chlorhexidine diacetate	13177008	2
Caelo	Chlorhexidine diacetate	13177021	8
Caelo	Chlorhexidine diacetate	1072015C	1
Caelo	Chlorhexidine diacetate	1317021	1
Caelo	Chlorhexidine diacetate	13177021	6
Caelo	Chlorhexidine diacetate	15261005	6
Caelo	Chlorhexidine diacetate	13177001	4
Caelo	Chlorhexidine diacetate	12142734	15
Caelo	Chlorhexidine diacetate	13177022	1
Caelo	Chlorhexidine diacetate	13177008	3
Caelo	Chlorhexidine diacetate	12177503	2
Caelo	Chlorhexidine diacetate	13177024	1
Caelo	Chlorhexidine diacetate	13177024	4
Caelo	Chlorhexidine diacetate	13177007	3
Caelo	Chlorhexidine diacetate	13177007	17
Caelo	Chlorhexidine diacetate	15438917	12
Caelo	Chlorhexidine diacetate	13177002	1
Caelo	Chlorhexidine diacetate	10154417	1
Caelo	Chlorhexidine diacetate	18072404	2
Caelo	Chlorhexidine diacetate	17325808	5
Caelo	Chlorhexidine diacetate	17326809	1
Caelo	Chlorhexidine diacetate	2195195	1
Caelo	Chlorhexidine diacetate	19121804	1
Caelo	Chlorhexidine diacetate	17325809	11
Caelo	Chlorhexidine diacetate	2195172	1
Caelo	Chlorhexidine diacetate	17325803	9
Caelo	Chlorhexidine diacetate	17325801	3
Caelo	Chlorhexidine diacetate	17185810	14
Caelo	Chlorhexidine diacetate	17185812	10
Caelo	Chlorhexidine diacetate	17188504	1
Caelo	Chlorhexidine diacetate	15438903	10
Caelo	Chlorhexidine diacetate	17185811	13
Caelo	Chlorhexidine diacetate	17185802	3
Caelo	Chlorhexidine diacetate	17185808	10
Caelo	Chlorhexidine diacetate	17185804	21
Caelo	Chlorhexidine diacetate	17185807	11
Caelo	Chlorhexidine diacetate	17185803	5
Caelo	Chlorhexidine diacetate	15438914	7
Caelo	Chlorhexidine diacetate	1543909	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Chlorhexidine diacetate	15438911	14
Caelo	Chlorhexidine diacetate	15438902	8
Caelo	Chlorhexidine diacetate	15438909	4
Caelo	Chlorhexidine diacetate	10154430	2
Caelo	Chlorhexidine diacetate	15261014	1
Caelo	Chlorhexidine diacetate	160912CH	1
Caelo	Chlorhexidine diacetate	180130CA/17185804	1
Caelo	Chlorhexidine diacetate	12142727	10
Caelo	Chlorhexidine diacetate	1735809	1
Caelo	Chlorhexidine diacetate	12177504	1
Caelo	Chlorhexidine diacetate	12142704	6
Caelo	Chlorhexidine diacetate	9113903	1
Caelo	Chlorhexidine diacetate	10154407	1
Caelo	Chlorhexidine diacetate	9113904	2
Caelo	Chlorhexidine diacetate	12142704ok	1
Caelo	Chlorhexidine diacetate	12142715	3
Caelo	Chlorhexidine diacetate	12142705	3
Caelo	Chlorhexidine diacetate	12142725	6
Caelo	Chlorhexidine diacetate	91139389	1
Caelo	Chlorhexidine diacetate	12142720	12
Caelo	Chlorhexidine diacetate	10154421	1
Caelo	Chlorhexidine diacetate	12177502	1
Caelo	Chlorhexidine diacetate	11204406	2
Caelo	Chlorhexidine diacetate	123456	1
Caelo	Chlorhexidine diacetate	12142709	1
Caelo	Chlorhexidine diacetate	12142726	10
Caelo	Chlorhexidine diacetate	11204401	3
Caelo	Chlorhexidine diacetate	1242720	1
Caelo	Chlorhexidine diacetate	11204407	1
Caelo	Chlorhexidine diacetate	12142702	2
Caelo	Chlorhexidine diacetate	10154429	1
Caelo	Chlorhexidine diacetate	12177510	9
Caelo	Chlorhexidine diacetate	10154424	1
Caelo	Chlorhexidine diacetate	11204405	1
Caelo	Chlorocresol	23415422	1
Caelo	Chlorocresol	4613M-01260	1
Caelo	Chlorocresol	15114107	1
Caelo	Chlorocresol	2041513	1
Caelo	Chlorocresol	10145239	3
Caelo	Chlorocresol	20101808	1
Caelo	Chlorocresol	12111013	2
Caelo	Chlorocresol	18027101	2
Caelo	Chlorocresol	7061304	1
Caelo	Chlorocresol	13275117	1
Caelo	Chlortetracycline hyd...	16B10-B09-347439	1
Caelo	Chlortetracycline hyd...	142008-03	1
Caelo	Chlortetracycline hyd...	FC15100204	5
Caelo	Chlortetracycline hyd...	14J08-B07-300249	1
Caelo	Chlortetracycline hyd...	FC13030134	1
Caelo	Chlortetracycline hyd...	160702Ü	1
Caelo	Chlortetracycline hyd...	141200803	2
Caelo	Cholesterol	11242103	1
Caelo	Cholesterol	12155605	1
Caelo	Cholesterol	15301803	3
Caelo	Cholesterol	12248403	1
Caelo	Cholesterol	17K16-B11-345814	1
Caelo	Cholesterol	16044415	1
Caelo	Cholesterol	27212	2
Caelo	Cholesterol	16G05-B03	1
Caelo	Cholesterol	13056903	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cholesterol	15407101	2
Caelo	Cholesterol	15407102	2
Caelo	Cholesterol	15J23-B03	1
Caelo	Cholesterol	15G31-B05-312376	1
Caelo	Cholesterol	15028303	1
Caelo	Cholesterol	14208602	1
Caelo	Cholesterol	14122501	2
Caelo	Cholesterol	13148704	1
Caelo	Cholesterol	14027201	1
Caelo	Cholesterol	15358801	5
Caelo	Ciclopirox olamine	14147402	2
Caelo	Ciclopirox olamine	12040906	1
Caelo	Ciclopirox olamine	13343	1
Caelo	Ciclopirox olamine	13343208	1
Caelo	Ciclopirox olamine	13343208	1
Caelo	Ciclopirox olamine	13343209	2
Caelo	Ciclopirox olamine	14147402	2
Caelo	Ciclopirox olamine	173074055727	1
Caelo	Ciclopirox olamine	15364409	1
Caelo	Ciclopirox olamine	15364404	1
Caelo	Ciclopirox olamine	15364407	4
Caelo	Ciclopirox olamine	15364408	4
Caelo	Ciclopirox olamine	17118811	2
Caelo	Ciclopirox olamine	14147405	1
Caelo	Ciclopirox olamine	15364412	6
Caelo	Ciclopirox olamine	15364405	3
Caelo	Ciclopirox olamine	13343201	1
Caelo	Ciclopirox olamine	15364416	2
Caelo	Ciclopirox olamine	17118801	2
Caelo	Ciclopirox olamine	16G26-B02-337355	1
Caelo	Ciclopirox olamine	17118802	3
Caelo	Ciclopirox olamine	17118807	2
Caelo	Ciclopirox olamine	17118808	5
Caelo	Ciclopirox olamine	15304407	1
Caelo	Ciclopirox olamine	14378802	1
Caelo	Ciclopirox olamine	15364402	2
Caelo	Ciclopirox olamine	17307408	1
Caelo	Ciclopirox olamine	14378810	4
Caelo	Ciclopirox olamine	14378805	2
Caelo	Ciclopirox olamine	14378806	5
Caelo	Ciclopirox olamine	14378803	1
Caelo	Ciclopirox olamine	18210201	2
Caelo	Ciclopirox olamine	18210202	1
Caelo	Ciclopirox olamine	14378801	2
Caelo	Ciclopirox olamine	15018803	1
Caelo	Ciclopirox olamine	17307402	2
Caelo	Ciclopirox olamine	13343203	1
Caelo	Ciclopirox olamine	15018802	2
Caelo	Ciclopirox olamine	14275801	2
Caelo	Ciclopirox olamine	14275802	2
Caelo	Ciclopirox olamine	13343209	1
Caelo	Ciclopirox olamine	17307405	2
Caelo	Ciclosporin	12165602	3
Caelo	Ciclosporin	12317703	1
Caelo	Ciclosporin	13148604	2
Caelo	Ciclosporin	13148606	3
Caelo	Ciclosporin	12317702	2
Caelo	Ciclosporin	13031303	1
Caelo	Ciclosporin	14199104	3
Caelo	Ciclosporin	12317707	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ciclosporin	13346808	1
Caelo	Ciclosporin	13346803	1
Caelo	Ciclosporin	14A15-No1	1
Caelo	Ciclosporin	13D26-N07	3
Caelo	Ciclosporin	13148606	1
Caelo	Ciclosporin	13346803	2
Caelo	Ciclosporin	19011301	1
Caelo	Ciclosporin	12165603	2
Caelo	Ciclosporin	17158309	1
Caelo	Ciclosporin	17332702	3
Caelo	Ciclosporin	12317707	3
Caelo	Ciclosporin	15278910	1
Caelo	Ciclosporin	19055804	1
Caelo	Ciclosporin	15278913	1
Caelo	Ciclosporin	15049703	3
Caelo	Ciclosporin	17158303	4
Caelo	Ciclosporin	18136108	2
Caelo	Ciclosporin	13148606	2
Caelo	Ciclosporin	17332703	1
Caelo	Ciclosporin	18136101	1
Caelo	Ciclosporin	15278902	2
Caelo	Ciclosporin	15278903	2
Caelo	Ciclosporin	29689114C22B02	1
Caelo	Ciclosporin	15278901	3
Caelo	Ciclosporin	17158306	1
Caelo	Ciclosporin	15278904	3
Caelo	Cinchona bark, cut	14009305	1
Caelo	Cinchona bark, cut	15301904	2
Caelo	Cinchona bark, cut	10381906	1
Caelo	Cinnamon rods (Ceylon...	12148904	1
Caelo	Cinnamon, cut	15151905	1
Caelo	Cinnamon, cut	12309908	1
Caelo	Cinnamon, cut	25111307	1
Caelo	Cinnamon, cut	13304401	2
Caelo	Cinnamon, cut	12309902	3
Caelo	Cinnamon, cut	15151902	1
Caelo	Cinnamon, cut	11357912	2
Caelo	Cinnamon, cut	13304409	1
Caelo	Cinnamon, cut	13304415	1
Caelo	Cinnamon, cut	13304402	2
Caelo	Cinnamon, cut	13304407	1
Caelo	Cinnamon, cut	16061702	1
Caelo	Ciprofloxacin hydroch...	16B17-B05-324218	1
Caelo	Citric acid	24021502	1
Caelo	Citric acid	13243005	1
Caelo	Citric acid	14186302	12
Caelo	Citric acid	1175409	2
Caelo	Citric acid	13396502	5
Caelo	Citric acid	13243216	12
Caelo	Citric acid	18091409	1
Caelo	Citric acid	13396505	3
Caelo	Citric acid	15091401	1
Caelo	Citric acid	13396504	3
Caelo	Citric acid	13243225	4
Caelo	Citric acid	13031502	1
Caelo	Citric acid	14041506	1
Caelo	Citric acid	14186301	26
Caelo	Citric acid	15063205	3
Caelo	Citric acid	15063211	4
Caelo	Citric acid	15087902	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Citric acid	15214301	1
Caelo	Citric acid	15063206	17
Caelo	Citric acid	15063202	1
Caelo	Citric acid	1044905	1
Caelo	Citric acid	15063003	23
Caelo	Citric acid	15044904	10
Caelo	Citric acid	14G04-B09-293876	8
Caelo	Citric acid	15063006	4
Caelo	Citric acid	15124301	8
Caelo	Citric acid	14071501	1
Caelo	Citric acid	15087510	1
Caelo	Citric acid	1072015B	1
Caelo	Citric acid	15044905	8
Caelo	Citric acid	15044902	4
Caelo	Citric acid	13243222	1
Caelo	Citric acid	14E06-B07	1
Caelo	Citric acid	133243218	1
Caelo	Citric acid	13061502	1
Caelo	Citric acid	15044901	6
Caelo	Citric acid	14186314	2
Caelo	Citric acid	14186313	9
Caelo	Citric acid	15044903	8
Caelo	Citric acid	13396506	1
Caelo	Citric acid	12353613	1
Caelo	Citric acid	11380213	1
Caelo	Citric acid	13243218	8
Caelo	Citric acid	12175708	3
Caelo	Citric acid	13G24-N13	10
Caelo	Citric acid	27091303	1
Caelo	Citric acid	13050623-3376319052	1
Caelo	Citric acid	10182420	1
Caelo	Citric acid	11380217	2
Caelo	Citric acid	13175405	5
Caelo	Citric acid	13050621	3
Caelo	Citric acid	12156822	2
Caelo	Citric acid	13175402	3
Caelo	Citric acid	13050622	7
Caelo	Citric acid	13175404	5
Caelo	Citric acid	17071301	1
Caelo	Citric acid	11C24-N03	1
Caelo	Citric acid	13050623	2
Caelo	Citric acid	13050610	6
Caelo	Citric acid	13050617	12
Caelo	Citric acid	2071303	1
Caelo	Citric acid	13050614	7
Caelo	Citric acid	1305601	1
Caelo	Citric acid	15063203	17
Caelo	Citric acid	12156830	1
Caelo	Citric acid	13050606	2
Caelo	Citric acid	12043119	1
Caelo	Citric acid	13050611	6
Caelo	Citric acid	1172673	11
Caelo	Citric acid	13K18-B34-288371TEST	1
Caelo	Citric acid	13175411	10
Caelo	Citric acid	13157404	1
Caelo	Citric acid	12175716	2
Caelo	Citric acid	13243224	2
Caelo	Citric acid	13243217	7
Caelo	Citric acid	160614	1
Caelo	Citric acid	13243221	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Citric acid	13243220	2
Caelo	Citric acid	13243210	9
Caelo	Citric acid	13243219	4
Caelo	Citric acid	13243208	2
Caelo	Citric acid	13243214	17
Caelo	Citric acid	28051401	1
Caelo	Citric acid	12175718	1
Caelo	Citric acid	13175409	21
Caelo	Citric acid	13175426	4
Caelo	Citric acid	13050604; 15R05	1
Caelo	Citric acid	13243202	11
Caelo	Citric acid	13243201	21
Caelo	Citric acid	13D29-N02	1
Caelo	Citric acid	28021405	1
Caelo	Citric acid	13175414	2
Caelo	Citric acid	13175406	3
Caelo	Citric acid	L130404202	1
Caelo	Citric acid	13050624	1
Caelo	Citric acid	12175711	2
Caelo	Citric acid	13175403	14
Caelo	Citric acid	26111301	1
Caelo	Citric acid	15087903	15
Caelo	Citric acid	16038713	1
Caelo	Citric acid	16004701	18
Caelo	Citric acid	17329504	1
Caelo	Citric acid	17162604	3
Caelo	Citric acid	1716262601	1
Caelo	Citric acid	18025702	30
Caelo	Citric acid	18025704	11
Caelo	Citric acid	5041802	1
Caelo	Citric acid	17329502	9
Caelo	Citric acid	162086055525	1
Caelo	Citric acid	18025801	14
Caelo	Citric acid	17162601	23
Caelo	Citric acid	13011801	1
Caelo	Citric acid	17162602	19
Caelo	Citric acid	16215409	4
Caelo	Citric acid	1716271	1
Caelo	Citric acid	16215405	4
Caelo	Citric acid	1620610	1
Caelo	Citric acid	17162603	4
Caelo	Citric acid	17162702	4
Caelo	Citric acid	18111704	1
Caelo	Citric acid	16038617	2
Caelo	Citric acid	16215408	3
Caelo	Citric acid	2617Q-01024	1
Caelo	Citric acid	16206813	1
Caelo	Citric acid	3091701	1
Caelo	Citric acid	1150817/16208614	1
Caelo	Citric acid	17162701	18
Caelo	Citric acid	16208615	3
Caelo	Citric acid	17329503	2
Caelo	Citric acid	16038712	4
Caelo	Citric acid	26051803	1
Caelo	Citric acid	18025701	20
Caelo	Citric acid	12156824	6
Caelo	Citric acid	18312104	1
Caelo	Citric acid	18306909	2
Caelo	Citric acid	18306905	4
Caelo	Citric acid	18306903	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Citric acid	18306906	4
Caelo	Citric acid	162154	1
Caelo	Citric acid	18097611	1
Caelo	Citric acid	18097604	1
Caelo	Citric acid	18306902	4
Caelo	Citric acid	17E15-B04-339371	1
Caelo	Citric acid	18097606	4
Caelo	Citric acid	18129206	1
Caelo	Citric acid	18129205	6
Caelo	Citric acid	18F25-B02-351453	1
Caelo	Citric acid	18097607	3
Caelo	Citric acid	18097609	9
Caelo	Citric acid	18F25-B02	2
Caelo	Citric acid	18129204	1
Caelo	Citric acid	17329512	1
Caelo	Citric acid	18097602	3
Caelo	Citric acid	17329513	2
Caelo	Citric acid	18097601	9
Caelo	Citric acid	18097603	10
Caelo	Citric acid	18025802	4
Caelo	Citric acid	180257045617	1
Caelo	Citric acid	17329511	3
Caelo	Citric acid	17161601	1
Caelo	Citric acid	16208612	17
Caelo	Citric acid	20170728/Ma1	1
Caelo	Citric acid	26071701	1
Caelo	Citric acid	16038608	11
Caelo	Citric acid	15063217	8
Caelo	Citric acid	15088008	3
Caelo	Citric acid	15088017	2
Caelo	Citric acid	150632	1
Caelo	Citric acid	16038707	10
Caelo	Citric acid	16038701	4
Caelo	Citric acid	1230616	1
Caelo	Citric acid	16004706	17
Caelo	Citric acid	13243229	1
Caelo	Citric acid	16038704	10
Caelo	Citric acid	1015U-01024	1
Caelo	Citric acid	15063003-20160615Ma1	1
Caelo	Citric acid	16038603	5
Caelo	Citric acid	16004705	3
Caelo	Citric acid	14K05-B04-302815	6
Caelo	Citric acid	16038604	13
Caelo	Citric acid	16004702	24
Caelo	Citric acid	15063213	25
Caelo	Citric acid	15088016	4
Caelo	Citric acid	16038602	3
Caelo	Citric acid	15063214	6
Caelo	Citric acid	15063218	10
Caelo	Citric acid	14G04B09293877	1
Caelo	Citric acid	1396505	1
Caelo	Citric acid	15,63218	1
Caelo	Citric acid	15124202	1
Caelo	Citric acid	21531	1
Caelo	Citric acid	160607AC	1
Caelo	Citric acid	14G04-B09	1
Caelo	Citric acid	16038613	9
Caelo	Citric acid	16215403	1
Caelo	Citric acid	6071713	1
Caelo	Citric acid	16208614	28

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Citric acid	16208613	8
Caelo	Citric acid	2317M-01024	1
Caelo	Citric acid	170607C/16208605	1
Caelo	Citric acid	16038614	4
Caelo	Citric acid	16208610	8
Caelo	Citric acid	16208605	12
Caelo	Citric acid	16208606	24
Caelo	Citric acid	16208607	26
Caelo	Citric acid	1917A-01024	1
Caelo	Citric acid	16038607	3
Caelo	Citric acid	15063204	9
Caelo	Citric acid	5315I-01024	1
Caelo	Citric acid	15185320	1
Caelo	Citric acid	16038622	5
Caelo	Citric acid	1501-B06-324348	1
Caelo	Citric acid	1516Q-01024	1
Caelo	Citric acid	16038709	5
Caelo	Citric acid	16038708	9
Caelo	Citric acid	8031702	1
Caelo	Citric acid	16038621	13
Caelo	Citric acid	xxxx	1
Caelo	Citric acid	xxxxxx	1
Caelo	Citric acid	15124303	1
Caelo	Citric acid	7808	1
Caelo	Citric acid	16038710	3
Caelo	Citric acid	2051703	1
Caelo	Citric acid	13050613	17
Caelo	Citric acid	13G24N13	2
Caelo	Citric acid	10049716	1
Caelo	Citric acid		4
Caelo	Citric acid	12156827	1
Caelo	Citric acid	11380215	3
Caelo	Citric acid	11159821	1
Caelo	Citric acid	9158408	1
Caelo	Citric acid	713I-01024	1
Caelo	Citric acid	12043110	5
Caelo	Citric acid	12043124	3
Caelo	Citric acid	12156803	8
Caelo	Citric acid	11066806	1
Caelo	Citric acid	11159801	2
Caelo	Citric acid	121568704	1
Caelo	Citric acid	20060912-3	1
Caelo	Citric acid	20151112-1	1
Caelo	Citric acid	1380208	1
Caelo	Citric acid	12043113	3
Caelo	Citric acid	11178303	5
Caelo	Citric acid	12156804	3
Caelo	Citric acid	12043120	5
Caelo	Citric acid	11K11-N02	3
Caelo	Citric acid	11189817	1
Caelo	Citric acid	11380204	5
Caelo	Citric acid	12156801	2
Caelo	Citric acid	9199105	1
Caelo	Citric acid	12156812	7
Caelo	Citric acid	19041305	1
Caelo	Citric acid	12043118	3
Caelo	Citric acid	11159824	1
Caelo	Citric acid	11380208	1
Caelo	Citric acid	A-040513-1	1
Caelo	Citric acid	16038606	19

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Citric acid	19041304	1
Caelo	Citric acid	13050607	10
Caelo	Citric acid	11159817	1
Caelo	Citric acid	13050601	10
Caelo	Citric acid	12156823	5
Caelo	Citric acid	12156818	1
Caelo	Citric acid	12043123	7
Caelo	Citric acid	12043105	2
Caelo	Citric acid	12043112	3
Caelo	Citric acid	13050604	15
Caelo	Citric acid monohydrate	14311807	2
Caelo	Citric acid monohydrate	14311808	4
Caelo	Citric acid monohydrate	11208132	3
Caelo	Citric acid monohydrate	sjdfksjfdkl	1
Caelo	Citric acid monohydrate	1215911	1
Caelo	Citric acid monohydrate	1215907	1
Caelo	Citric acid monohydrate	12125907	6
Caelo	Citric acid monohydrate	14038403	7
Caelo	Citric acid monohydrate	12125938	1
Caelo	Citric acid monohydrate	13069133	5
Caelo	Citric acid monohydrate	13069104	6
Caelo	Citric acid monohydrate	13069129	4
Caelo	Citric acid monohydrate	12125937	4
Caelo	Citric acid monohydrate	1212937	1
Caelo	Citric acid monohydrate	13069124	4
Caelo	Citric acid monohydrate	14038401	8
Caelo	Citric acid monohydrate	12125904	2
Caelo	Citric acid monohydrate	13069106	3
Caelo	Citric acid monohydrate	13069130	1
Caelo	Citric acid monohydrate	13069105	2
Caelo	Citric acid monohydrate	13069110	4
Caelo	Citric acid monohydrate	13069111	4
Caelo	Citric acid monohydrate	13069113	2
Caelo	Citric acid monohydrate	13069128	3
Caelo	Citric acid monohydrate	13069125	2
Caelo	Citric acid monohydrate	11208133	4
Caelo	Citric acid monohydrate	16186501	2
Caelo	Citric acid monohydrate	17085607	4
Caelo	Citric acid monohydrate	17329804	4
Caelo	Citric acid monohydrate	17085603	3
Caelo	Citric acid monohydrate	17085602	8
Caelo	Citric acid monohydrate	16186508	5
Caelo	Citric acid monohydrate	17329807	1
Caelo	Citric acid monohydrate	14071701	1
Caelo	Citric acid monohydrate	16186502	5
Caelo	Citric acid monohydrate	17329801	1
Caelo	Citric acid monohydrate	17329808	1
Caelo	Citric acid monohydrate	14311812	4
Caelo	Citric acid monohydrate	15185320	4
Caelo	Citric acid monohydrate	15185321	3
Caelo	Citric acid monohydrate	15185316	5
Caelo	Citric acid monohydrate	19088402	1
Caelo	Citric acid monohydrate	18130503	2
Caelo	Citric acid monohydrate	15185317	1
Caelo	Citric acid monohydrate	15087502	2
Caelo	Citric acid monohydrate	15087503	2
Caelo	Citric acid monohydrate	15087509	8
Caelo	Citric acid monohydrate	13069103	1
Caelo	Citric acid monohydrate	15087510	2
Caelo	Citric acid monohydrate	14311813	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Citric acid monohydrate	15185314	5
Caelo	Citric acid monohydrate	15087512	1
Caelo	Citric acid monohydrate	2515Q-02959	1
Caelo	Citric acid monohydrate	15185302	11
Caelo	Citric acid monohydrate	15185301	2
Caelo	Citric acid monohydrate	15185309	2
Caelo	Citric acid monohydrate	15185304	10
Caelo	Citric acid monohydrate	161865001	1
Caelo	Clindamycin hydrochlo...	1811003-01	4
Caelo	Clindamycin hydrochlo...		1
Caelo	Clindamycin hydrochlo...	16H03-B03-334345	1
Caelo	Clindamycin hydrochlo...	17001301	1
Caelo	Clindamycin hydrochlo...	16H03-B02	2
Caelo	Clindamycin hydrochlo...	17209303	2
Caelo	Clindamycin hydrochlo...	17209302	3
Caelo	Clindamycin hydrochlo...	18J15-B04-355581	1
Caelo	Clindamycin hydrochlo...	160712C	1
Caelo	Clindamycin hydrochlo...	10809406	1
Caelo	Clindamycin hydrochlo...	18A18-B08-346106	3
Caelo	Clindamycin hydrochlo...	18009401	1
Caelo	Clindamycin hydrochlo...	15392822	2
Caelo	Clindamycin hydrochlo...	18009413	1
Caelo	Clindamycin hydrochlo...	P-003-WS10100301	1
Caelo	Clindamycin hydrochlo...	18009405	2
Caelo	Clindamycin hydrochlo...	151003801	2
Caelo	Clindamycin hydrochlo...	15392802	1
Caelo	Clindamycin hydrochlo...	15392815	4
Caelo	Clindamycin hydrochlo...	15L22-B03-318788	1
Caelo	Clindamycin hydrochlo...	15392807	1
Caelo	Clindamycin hydrochlo...	15392816	1
Caelo	Clindamycin hydrochlo...	15392804	3
Caelo	Clindamycin hydrochlo...	15392811	2
Caelo	Clindamycin hydrochlo...	1511004-02	5
Caelo	Clindamycin hydrochlo...	15392805	1
Caelo	Clindamycin hydrochlo...	15392817	2
Caelo	Clindamycin hydrochlo...	15392814	1
Caelo	Clindamycin hydrochlo...	15392818	2
Caelo	Clindamycin hydrochlo...	15189512	1
Caelo	Clindamycin hydrochlo...	14K04B05309851	1
Caelo	Clindamycin hydrochlo...	15392819	1
Caelo	Clindamycin-2-dihydro...	14257305	1
Caelo	Clindamycin-2-dihydro...	17344101	1
Caelo	Clioquinol	14359117	2
Caelo	Clioquinol	18011501	1
Caelo	Clioquinol	1214I02562	1
Caelo	Clioquinol	13114413	10
Caelo	Clioquinol	18011516	10
Caelo	Clioquinol	1509038-OTC-Pharma	1
Caelo	Clioquinol	16307215	7
Caelo	Clioquinol	180115117	1
Caelo	Clioquinol	18011523	2
Caelo	Clioquinol	11(18011523)	1
Caelo	Clioquinol	14359118	14
Caelo	Clioquinol	14359113	2
Caelo	Clioquinol	18278206	1
Caelo	Clioquinol	14359120	8
Caelo	Clioquinol	170926CL/16307208	1
Caelo	Clioquinol	18011502	25
Caelo	Clioquinol	11304102	7
Caelo	Clioquinol	3513I-02562	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clioquinol	11304119	4
Caelo	Clioquinol	11304107	6
Caelo	Clioquinol	11304120	14
Caelo	Clioquinol	11304112	7
Caelo	Clioquinol	13114404	2
Caelo	Clioquinol	11304103	1
Caelo	Clioquinol	11304117	10
Caelo	Clioquinol	150903803	1
Caelo	Clioquinol	13A29-N07	8
Caelo	Clioquinol	13114408	5
Caelo	Clioquinol	11304109	1
Caelo	Clioquinol	11304113	4
Caelo	Clioquinol	10341608	1
Caelo	Clioquinol	1130312	1
Caelo	Clioquinol	14/-B03-309181	1
Caelo	Clioquinol	130114-02	1
Caelo	Clioquinol	4613E2562	1
Caelo	Clioquinol	144359111	1
Caelo	Clioquinol	11304124	15
Caelo	Clioquinol	5213A2562	1
Caelo	Clioquinol	11304116	3
Caelo	Clioquinol	13114406	1
Caelo	Clioquinol	5031302	1
Caelo	Clioquinol	14359112	13
Caelo	Clioquinol	103415	2
Caelo	Clioquinol	16307210	2
Caelo	Clioquinol	16307207	16
Caelo	Clioquinol	13114422	8
Caelo	Clioquinol	5089	1
Caelo	Clioquinol	16307208	9
Caelo	Clioquinol	13114424	11
Caelo	Clioquinol	1201034-01	7
Caelo	Clioquinol	14359111	13
Caelo	Clioquinol	16118-B04-329141	1
Caelo	Clioquinol	10341601	2
Caelo	Clioquinol	16307203	2
Caelo	Clioquinol	17204603	9
Caelo	Clioquinol	16307213	1
Caelo	Clioquinol	13114423	10
Caelo	Clioquinol	16307209	1
Caelo	Clioquinol	16307216	8
Caelo	Clioquinol	10341605	4
Caelo	Clioquinol	16021505	1
Caelo	Clioquinol	1405080031	1
Caelo	Clioquinol	13114421	1
Caelo	Clioquinol	13114427	2
Caelo	Clioquinol	163072	1
Caelo	Clioquinol	16307214	10
Caelo	Clioquinol	16307211	6
Caelo	Clioquinol	14359102	10
Caelo	Clioquinol	1405080036	1
Caelo	Clioquinol	11304105	3
Caelo	Clioquinol	13114417	10
Caelo	Clioquinol	16122020-01	1
Caelo	Clioquinol	4214Q-02562	1
Caelo	Clioquinol	14/11-BG03-309181	1
Caelo	Clioquinol	14359119	3
Caelo	Clioquinol	3041409	1
Caelo	Clioquinol	13114410	2
Caelo	Clioquinol	17204601	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clioquinol	14353120	1
Caelo	Clioquinol	10341510	1
Caelo	Clioquinol	160916CL	1
Caelo	Clioquinol	3713A2562	1
Caelo	Clioquinol	13114415	3
Caelo	Clioquinol	14359116	1
Caelo	Clioquinol	180127C/16307215	1
Caelo	Clioquinol	16307204	1
Caelo	Clioquinol	16307206	19
Caelo	Clioquinol	16307205	8
Caelo	Clobetasol propionate	1524003	1
Caelo	Clobetasol propionate	14338006	12
Caelo	Clobetasol propionate	CL0/1191217	1
Caelo	Clobetasol propionate	18082702	2
Caelo	Clobetasol propionate	14338008	4
Caelo	Clobetasol propionate	15234002	3
Caelo	Clobetasol propionate	8082703	1
Caelo	Clobetasol propionate	15234004	23
Caelo	Clobetasol propionate	18/12-B02-35480	1
Caelo	Clobetasol propionate	17325605	10
Caelo	Clobetasol propionate	1105044-03	1
Caelo	Clobetasol propionate	18082703	10
Caelo	Clobetasol propionate	13K08-B30-287726	3
Caelo	Clobetasol propionate	15234001	12
Caelo	Clobetasol propionate	2x17325602	1
Caelo	Clobetasol propionate	17325606	3
Caelo	Clobetasol propionate	13134703	4
Caelo	Clobetasol propionate	13269507	4
Caelo	Clobetasol propionate	15A14-B02	1
Caelo	Clobetasol propionate	14C03-B07-293560	3
Caelo	Clobetasol propionate	14338002	1
Caelo	Clobetasol propionate	15439504	25
Caelo	Clobetasol propionate	15412910	1
Caelo	Clobetasol propionate	14338007	6
Caelo	Clobetasol propionate	15234003	20
Caelo	Clobetasol propionate	17096701	10
Caelo	Clobetasol propionate	13269502	3
Caelo	Clobetasol propionate	16C18-B02-325573	1
Caelo	Clobetasol propionate	14338001	9
Caelo	Clobetasol propionate	15439502	7
Caelo	Clobetasol propionate	13K08-B30-287728	6
Caelo	Clobetasol propionate	17096704	1
Caelo	Clobetasol propionate	16106718	1
Caelo	Clobetasol propionate	11101311	1
Caelo	Clobetasol propionate	18082718	1
Caelo	Clobetasol propionate	14188506	5
Caelo	Clobetasol propionate	14g30bo7302114	1
Caelo	Clobetasol propionate	1438007	1
Caelo	Clobetasol propionate	16106719	4
Caelo	Clobetasol propionate	14338005	3
Caelo	Clobetasol propionate	18082714	1
Caelo	Clobetasol propionate	apl/294/b-14	1
Caelo	Clobetasol propionate	14188503	15
Caelo	Clobetasol propionate	2014E-02835	1
Caelo	Clobetasol propionate	16353208	1
Caelo	Clobetasol propionate	17325603	14
Caelo	Clobetasol propionate	14188502	7
Caelo	Clobetasol propionate	13269509	1
Caelo	Clobetasol propionate	13134702	7
Caelo	Clobetasol propionate	14102416	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clobetasol propionate	18082708	8
Caelo	Clobetasol propionate	14188504	7
Caelo	Clobetasol propionate	1616715	1
Caelo	Clobetasol propionate	17090708	1
Caelo	Clobetasol propionate	1369504	1
Caelo	Clobetasol propionate	14188501	3
Caelo	Clobetasol propionate	13269508	13
Caelo	Clobetasol propionate	16106712	3
Caelo	Clobetasol propionate	1704020	2
Caelo	Clobetasol propionate	1606022-025349	1
Caelo	Clobetasol propionate	12342312	1
Caelo	Clobetasol propionate	17096707	4
Caelo	Clobetasol propionate	1704020-01;180522CP	1
Caelo	Clobetasol propionate	13134701	15
Caelo	Clobetasol propionate	17096702	6
Caelo	Clobetasol propionate	1735603	1
Caelo	Clobetasol propionate	13134704	4
Caelo	Clobetasol propionate	14G30_B07_302114	1
Caelo	Clobetasol propionate	15031401	1
Caelo	Clobetasol propionate	17325602	13
Caelo	Clobetasol propionate	13134705	5
Caelo	Clobetasol propionate	13269501	1
Caelo	Clobetasol propionate	17096709	6
Caelo	Clobetasol propionate	152340004	1
Caelo	Clobetasol propionate	12342307	2
Caelo	Clobetasol propionate	12107609	7
Caelo	Clobetasol propionate	1408018-03	32
Caelo	Clobetasol propionate	16106702	12
Caelo	Clobetasol propionate	11300707	2
Caelo	Clobetasol propionate	170404CP1606022-02	1
Caelo	Clobetasol propionate	16106703	18
Caelo	Clobetasol propionate	181016709	1
Caelo	Clobetasol propionate	15j15-b04-317698	1
Caelo	Clobetasol propionate	15253032	1
Caelo	Clobetasol propionate	12107604	1
Caelo	Clobetasol propionate	12107602	3
Caelo	Clobetasol propionate	20121012-4	1
Caelo	Clobetasol propionate	15439507	2
Caelo	Clobetasol propionate	123645	1
Caelo	Clobetasol propionate	16E04-B12	2
Caelo	Clobetasol propionate	23111202	1
Caelo	Clobetasol propionate	31071302	1
Caelo	Clobetasol propionate	12107601	3
Caelo	Clobetasol propionate	11300702	4
Caelo	Clobetasol propionate	15439501	7
Caelo	Clobetasol propionate	11101315	1
Caelo	Clobetasol propionate	16106701	3
Caelo	Clobetasol propionate	16106709	13
Caelo	Clobetasol propionate	1616703	1
Caelo	Clobetasol propionate	16106710	4
Caelo	Clobetasol propionate	16H22-B07-328012	12
Caelo	Clobetasol propionate	18H22-B017-328012	1
Caelo	Clobetasol propionate	16106704	23
Caelo	Clobetasol propionate	154390504	1
Caelo	Clobetasol propionate	3012Q02835	1
Caelo	Clobetasol propionate	12107607	6
Caelo	Clobetasol propionate	CL0024/0915	1
Caelo	Clobetasol propionate	16348001	1
Caelo	Clobetasol propionate	12342303	2
Caelo	Clobetasol propionate	16106713	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clobetasol propionate	15234009	5
Caelo	Clobetasol propionate	12342306	9
Caelo	Clobetasol propionate	16106715	9
Caelo	Clobetasol propionate	15234008	1
Caelo	Clobetasol propionate	15J-15-B04-316616	1
Caelo	Clobetasol propionate	12342309	13
Caelo	Clobetasol propionate	15439503	22
Caelo	Clobetasol propionate	16106711	17
Caelo	Clobetasol propionate	160115CP	1
Caelo	Clobetasol propionate	160829CP	1
Caelo	Clobetasol propionate	1610670	1
Caelo	Clotrimazole	17360204	8
Caelo	Clotrimazole	16234504	11
Caelo	Clotrimazole	132243218	1
Caelo	Clotrimazole	14144905	18
Caelo	Clotrimazole	18147203	14
Caelo	Clotrimazole	18147206	1
Caelo	Clotrimazole	18147103	11
Caelo	Clotrimazole	18147101	13
Caelo	Clotrimazole	181471035763	1
Caelo	Clotrimazole	18A18-B05-349662	1
Caelo	Clotrimazole	18147210	8
Caelo	Clotrimazole	700135-0006	1
Caelo	Clotrimazole	18147211	5
Caelo	Clotrimazole	17260204	1
Caelo	Clotrimazole	18081408	1
Caelo	Clotrimazole	19016108	1
Caelo	Clotrimazole	12B27-N05	4
Caelo	Clotrimazole	12320806	2
Caelo	Clotrimazole	12320809	4
Caelo	Clotrimazole	12320808	11
Caelo	Clotrimazole	1259913	1
Caelo	Clotrimazole	18147208	1
Caelo	Clotrimazole	1234567	1
Caelo	Clotrimazole	1711604	2
Caelo	Clotrimazole	171113608	1
Caelo	Clotrimazole	16J07-B02_330236	1
Caelo	Clotrimazole	14144912	3
Caelo	Clotrimazole	B: 1024305	1
Caelo	Clotrimazole	11-K21-N19	1
Caelo	Clotrimazole	14144911	12
Caelo	Clotrimazole	10247328	3
Caelo	Clotrimazole	1105014-02	2
Caelo	Clotrimazole	11/22-N10	1
Caelo	Clotrimazole	17111617	2
Caelo	Clotrimazole	16B11-B04-322393	1
Caelo	Clotrimazole	11305506	2
Caelo	Clotrimazole	12159902	12
Caelo	Clotrimazole	20111115	2
Caelo	Clotrimazole	12159904	2
Caelo	Clotrimazole	11305511	1
Caelo	Clotrimazole	12159911	2
Caelo	Clotrimazole	17360202	14
Caelo	Clotrimazole		5
Caelo	Clotrimazole	14337824NP	1
Caelo	Clotrimazole	17111603	3
Caelo	Clotrimazole	12159916	6
Caelo	Clotrimazole	12159913	9
Caelo	Clotrimazole	12320807	4
Caelo	Clotrimazole	12320813	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clotrimazole	17/25-B06-345048	1
Caelo	Clotrimazole	12320805	3
Caelo	Clotrimazole	12320814	11
Caelo	Clotrimazole	12320815	2
Caelo	Clotrimazole	13148007	1
Caelo	Clotrimazole	17111615	14
Caelo	Clotrimazole	1306036	2
Caelo	Clotrimazole	1711613	1
Caelo	Clotrimazole	13224204	2
Caelo	Clotrimazole	1711816	1
Caelo	Clotrimazole	1711606	1
Caelo	Clotrimazole	1711603	2
Caelo	Clotrimazole	17111605	11
Caelo	Clotrimazole	1309034-C1	1
Caelo	Clotrimazole	14144913	3
Caelo	Clotrimazole	14337801	14
Caelo	Clotrimazole	13224215	4
Caelo	Clotrimazole	156546	1
Caelo	Clotrimazole	17111601	7
Caelo	Clotrimazole	132234216	1
Caelo	Clotrimazole	13324204	1
Caelo	Clotrimazole	17111616	11
Caelo	Clotrimazole	14144917	5
Caelo	Clotrimazole	16J07-B02	2
Caelo	Clotrimazole	14144918	12
Caelo	Clotrimazole	13224220	6
Caelo	Clotrimazole	13224217	8
Caelo	Clotrimazole	415050-2-P	1
Caelo	Clotrimazole	Fagron	1
Caelo	Clotrimazole	17111613	18
Caelo	Clotrimazole	160614	1
Caelo	Clotrimazole	5117E-01276	1
Caelo	Clotrimazole	13224216	10
Caelo	Clotrimazole	13184103	8
Caelo	Clotrimazole	26071305	1
Caelo	Clotrimazole	1882321	1
Caelo	Clotrimazole	11305507	1
Caelo	Clotrimazole	13184111	7
Caelo	Clotrimazole	13224201	26
Caelo	Clotrimazole	13184113	2
Caelo	Clotrimazole	311111	1
Caelo	Clotrimazole	1210030-02	1
Caelo	Clotrimazole	131841111	1
Caelo	Clotrimazole	13224219	7
Caelo	Clotrimazole	14c04-b30-291656	1
Caelo	Clotrimazole	14337802	7
Caelo	Clotrimazole	13224210	12
Caelo	Clotrimazole	PC85	1
Caelo	Clotrimazole	1322401	1
Caelo	Clotrimazole	17111606	13
Caelo	Clotrimazole	13224211	5
Caelo	Clotrimazole	418072-1-P	1
Caelo	Clotrimazole	11081201	1
Caelo	Clotrimazole	13224212	8
Caelo	Clotrimazole	162345225555	1
Caelo	Clotrimazole	13224221	2
Caelo	Clotrimazole	412091-2-P	1
Caelo	Clotrimazole	13224218	11
Caelo	Clotrimazole	14144914	9
Caelo	Clotrimazole	1711161	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clotrimazole	14337809	13
Caelo	Clotrimazole	15G14-B06-311057	1
Caelo	Clotrimazole	14337814	9
Caelo	Clotrimazole	16234503	12
Caelo	Clotrimazole	14C04-B31-299955	1
Caelo	Clotrimazole	15G14-B05-316483	4
Caelo	Clotrimazole	15395603	9
Caelo	Clotrimazole	14337828	7
Caelo	Clotrimazole	14337827	11
Caelo	Clotrimazole	14337832	1
Caelo	Clotrimazole	15395602	15
Caelo	Clotrimazole	16234520	14
Caelo	Clotrimazole	14337815	12
Caelo	Clotrimazole	14c04-b31-299954	2
Caelo	Clotrimazole	Val. Mes : 2815M01276	1
Caelo	Clotrimazole	14337819	7
Caelo	Clotrimazole	15295603	1
Caelo	Clotrimazole	14337830	2
Caelo	Clotrimazole	14C04B31-299954	1
Caelo	Clotrimazole	14337813	6
Caelo	Clotrimazole	15395604	13
Caelo	Clotrimazole	14337816	15
Caelo	Clotrimazole	16234523	1
Caelo	Clotrimazole	16234512	16
Caelo	Clotrimazole	14c04-B31	1
Caelo	Clotrimazole	17111604	6
Caelo	Clotrimazole	15G14-B05-311101	2
Caelo	Clotrimazole	153965604	1
Caelo	Clotrimazole	162345135319	1
Caelo	Clotrimazole	15395615	3
Caelo	Clotrimazole	1437816	1
Caelo	Clotrimazole	15395613	9
Caelo	Clotrimazole	16J07-B02-329385	4
Caelo	Clotrimazole	143337811	1
Caelo	Clotrimazole	16234506	1
Caelo	Clotrimazole	16B09B06311506	1
Caelo	Clotrimazole	15395623	2
Caelo	Clotrimazole	16234513	16
Caelo	Clotrimazole	16B11-B04-325221	1
Caelo	Clotrimazole	16234505	19
Caelo	Clotrimazole	16234502	4
Caelo	Clotrimazole	141104C	1
Caelo	Clotrimazole	15395618	16
Caelo	Clotrimazole	16234509	11
Caelo	Clotrimazole	16E20-B10	2
Caelo	Clotrimazole	14337818	12
Caelo	Clotrimazole	15395611	8
Caelo	Clotrimazole	15395614	9
Caelo	Clotrimazole	16234522	8
Caelo	Clotrimazole	1539613	1
Caelo	Clotrimazole	16E20-B10-324555	4
Caelo	Clotrimazole	16234517	6
Caelo	Clotrimazole	16234525	2
Caelo	Clotrimazole	17053186	1
Caelo	Clotrimazole	14337817	15
Caelo	Clotrimazole	1623455	1
Caelo	Clotrimazole	310316-5-461	1
Caelo	Clotrimazole	16051913	1
Caelo	Clotrimazole	14337811	5
Caelo	Clotrimazole	14337824	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Clotrimazole	14337829	16
Caelo	Clotrimazole	1170717	1
Caelo	Clotrimazole	14337812	7
Caelo	Clotrimazole	2217Q-01276	1
Caelo	Clotrimazole	17011601	1
Caelo	Clotrimazole	14337807	12
Caelo	Clotrimazole	15395605	11
Caelo	Clove	13371309	1
Caelo	Clove	14057602	1
Caelo	Clove	11233305	1
Caelo	Clove	12095806	2
Caelo	Clove	12095807	5
Caelo	Clove	12266505	1
Caelo	Clove	12095802	1
Caelo	Clove	5398872	1
Caelo	Clove, powdered (spice)	12355303	1
Caelo	Clove, powdered (spice)	14180001	1
Caelo	Cocoa bark, cut	12126907	1
Caelo	Cocoa bark, cut	19011303	1
Caelo	Cocoa bark, cut	12126905	2
Caelo	Cocoa bark, cut	1216907	1
Caelo	Cocoa bark, cut	12126910	1
Caelo	Cocoa bark, cut	12126901	1
Caelo	Cocoa butter	13292013	9
Caelo	Cocoa butter	13292015	1
Caelo	Cocoa butter	1392013	1
Caelo	Cocoa butter	18135102	2
Caelo	Cocoa butter	13292019	3
Caelo	Cocoa butter	14123611	6
Caelo	Cocoa butter	18135106	5
Caelo	Cocoa butter	13140409	1
Caelo	Cocoa butter	14123617	3
Caelo	Cocoa butter	16316509	5
Caelo	Cocoa butter	15246102	4
Caelo	Cocoa butter	14123609	2
Caelo	Cocoa butter	14045104	1
Caelo	Cocoa butter	515A-03176	1
Caelo	Cocoa butter	14348908	2
Caelo	Cocoa butter	14304703	1
Caelo	Cocoa butter	15088904	1
Caelo	Cocoa butter	14348906	8
Caelo	Cocoa butter	18135107	1
Caelo	Cocoa butter	14123620	5
Caelo	Cocoa butter	18135101	1
Caelo	Cocoa butter	15088903	1
Caelo	Cocoa butter	13292001	2
Caelo	Cocoa butter	17314901	2
Caelo	Cocoa butter	16049209	2
Caelo	Cocoa butter	13027207	2
Caelo	Cocoa butter	163116508	1
Caelo	Cocoa butter	13027221	4
Caelo	Cocoa butter	17055902	9
Caelo	Cocoa butter	13027224	5
Caelo	Cocoa butter	17055907	12
Caelo	Cocoa butter	17055905	1
Caelo	Cocoa butter	13027226	4
Caelo	Cocoa butter	13140401	1
Caelo	Cocoa butter	16062604	1
Caelo	Cocoa butter	12270209	2
Caelo	Cocoa butter	16316501	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Cocoa butter	170555901	1
Caelo	Cocoa butter	16049201	2
Caelo	Cocoa butter	17314902	17
Caelo	Cocoa butter	14348909	2
Caelo	Cocoa butter	13140404	3
Caelo	Cocoa butter	13027227	3
Caelo	Cocoa butter	13292006	4
Caelo	Cocoa butter	13292007	10
Caelo	Cocoa butter	13140405	3
Caelo	Cocoa butter	13292011	3
Caelo	Cocoa butter	13140407	4
Caelo	Cocoa butter	13140406	3
Caelo	Cocoa butter	13292012	3
Caelo	Cocoa butter	17055913	1
Caelo	Cocoa butter	16316503	6
Caelo	Cocoa butter	14123621	2
Caelo	Cocoa butter	12270202	4
Caelo	Cocoa butter	12270201	2
Caelo	Cocoa butter	11117816	2
Caelo	Cocoa butter	103276705	1
Caelo	Cocoa butter	12342003	7
Caelo	Cocoa butter	12040802	1
Caelo	Cocoa butter	12270203	4
Caelo	Cocoa butter	13027204	4
Caelo	Cocoa butter	13027202	4
Caelo	Cocoa butter	11034719	1
Caelo	Cocoa butter	13027203	1
Caelo	Cocoa butter	11034718	1
Caelo	Cocoa butter	12040807	2
Caelo	Cocoa butter	12040803	1
Caelo	Cocoa butter	18135112	1
Caelo	Colistin sulfate	1602021-01	1
Caelo	Coltsfoot leaf, cut	23041302	1
Caelo	Coltsfoot leaf, cut	Int.Ch:Now16081	1
Caelo	Comfrey root	116053	1
Caelo	Comfrey root	2711Q-02135	1
Caelo	Common broom herb, cut	13241601	1
Caelo	Common broom herb, cut	12358805	1
Caelo	Common broom herb, cut	4412Q-01680	1
Caelo	Common broom herb, cut	11132209	1
Caelo	Common broom herb, cut	12358804	3
Caelo	Common broom herb, cut	12358802	1
Caelo	Common wormwood herb,...	12016709	2
Caelo	Common wormwood herb,...	14360605	1
Caelo	Common wormwood herb,...	12016707	3
Caelo	Common wormwood herb,...	12382705	2
Caelo	Common wormwood herb,...	10355329	1
Caelo	Common wormwood herb,...	12382701	7
Caelo	Common wormwood herb,...	12382702	12
Caelo	Common wormwood herb,...	12016702	1
Caelo	Common wormwood herb,...	18031302	1
Caelo	Common wormwood herb,...	12016705	5
Caelo	Common wormwood herb,...	1081301	1
Caelo	Common wormwood herb,...	12382706	1
Caelo	Common wormwood herb,...	9355403	1
Caelo	Common wormwood herb,...	12016711	3
Caelo	Common wormwood herb,...	1238701	1
Caelo	Common wormwood herb,...	3591A110119	4
Caelo	Common wormwood herb,...	12016708	5
Caelo	Copper (H1032)	933901	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Copper (II) sulfate p...	12157107	1
Caelo	Copper (II) sulfate p...	3211M-01317	1
Caelo	Copper (II) sulfate p...	10L15-N06	1
Caelo	Copper (II) sulfate p...	2913U-01317	1
Caelo	Copper (II) sulfate p...	3411M-01317	1
Caelo	Coriander	12080906	1
Caelo	Coriander	14031316	1
Caelo	Coriander	12080905	1
Caelo	Coriander	13010306	2
Caelo	Coriander	4091305	1
Caelo	Coriander	16010209	1
Caelo	Coriander	9285707	1
Caelo	Coriander	27061305	1
Caelo	Coriander	13010302	2
Caelo	Coriander	13010304	1
Caelo	Coriander	16010201	1
Caelo	Coriander	18010209	1
Caelo	Coriander	12080901	2
Caelo	Corn starch	12238715	2
Caelo	Corn starch	13154904	1
Caelo	Corn starch	15119011	5
Caelo	Corn starch	13154914	1
Caelo	Corn starch	15310903	2
Caelo	Corn starch	11179307	2
Caelo	Corn starch	15310914	5
Caelo	Corn starch	15310909	2
Caelo	Corn starch	16138505	4
Caelo	Corn starch	16138507	2
Caelo	Corn starch	14111701	1
Caelo	Corn starch	17047404	6
Caelo	Corn starch	15310904	1
Caelo	Corn starch	10061146	1
Caelo	Corn starch	18274802	1
Caelo	Corn starch	12238705	1
Caelo	Cottonwood bud	70864257	1
Caelo	Cottonwood bud	12027504	1
Caelo	Cottonwood bud	12027506	1
Caelo	Couch grass rhizome, ...	10349705	1
Caelo	Couch grass rhizome, ...	3712I-02180	1
Caelo	Couch grass rhizome, ...	10349704	1
Caelo	Couch grass rhizome, ...	10049614	1
Caelo	Couch grass rhizome, ...	12186817	3
Caelo	Couch grass rhizome, ...	10049618	2
Caelo	Couch grass rhizome, ...	12186705	1
Caelo	Cowberry leaf, cut	11334302	1
Caelo	Cowberry leaf, cut	12105203	1
Caelo	Cowberry leaf, cut	12150205	1
Caelo	Cowberry leaf, cut	12105201	1
Caelo	Cyanocobalamin	9373503	2
Caelo	Cyanocobalamin	9324113	1
Caelo	Cyanocobalamin	9351704	1
Caelo	Cyanocobalamin	1604017-01	1
Caelo	Cyproterone acetate	17D11-B07-334581	3
Caelo	Cyproterone acetate	13j16-n04	1
Caelo	Daisy flower	12133703	1
Caelo	Daisy flower	11103203	2
Caelo	Daisy flower	12133704	1
Caelo	Daisy flower	1210a120215	2
Caelo	Daisy flower	12133705	1
Caelo	Daisy flower	11103201	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Daisy flower	11068308	1
Caelo	Damiana leaf, cut	12330402	7
Caelo	Damiana leaf, cut	11026607	2
Caelo	Damiana leaf, cut	12330401	5
Caelo	Damiana leaf, cut	12064103	3
Caelo	Damiana leaf, cut	12330403	2
Caelo	Damiana leaf, cut	12330407	2
Caelo	Dandelion herb with r...	17121305	1
Caelo	Dandelion herb with r...	11030314	4
Caelo	Dandelion herb with r...	13199202	2
Caelo	Dandelion herb with r...	Int.Ch:318090	1
Caelo	Dandelion herb with r...	12236002	8
Caelo	Dandelion herb with r...	13317903	1
Caelo	Dandelion herb with r...	2112122165	1
Caelo	Dandelion herb with r...	11030309	2
Caelo	Dandelion herb with r...	11030311	2
Caelo	Dandelion herb with r...	11030305	1
Caelo	Dandelion herb with r...	6611B100902	1
Caelo	Dandelion herb with r...	20111202	1
Caelo	Dandelion herb with r...	11371903	1
Caelo	Dandelion herb with r...	7341703	1
Caelo	Dandelion herb with r...	1103011	1
Caelo	Dandelion herb with r...	13108801	4
Caelo	Dandelion herb with r...	13018902	1
Caelo	Dandelion herb with r...	12236001	3
Caelo	Dandelion herb with r...	Ch.-B.:13019802	1
Caelo	Dandelion herb with r...	806212165	1
Caelo	Dandelion herb with r...	12236005	5
Caelo	Dandelion herb with r...	13019805	1
Caelo	Dandelion herb with r...	11030313	2
Caelo	Dandelion herb with r...	12236003	3
Caelo	Dandelion herb with r...	13108802	11
Caelo	Dandelion herb with r...	10121602	1
Caelo	Dandelion herb with r...	8111205	1
Caelo	Dandelion herb with r...	13019802	19
Caelo	Dandelion herb with r...	13019801	5
Caelo	Dandelion herb with r...	13019802-3395310101	1
Caelo	Dandelion herb withou...	1127102	1
Caelo	Dandelion herb withou...	9011403	1
Caelo	Dandelion herb withou...	107063	1
Caelo	Dandelion herb withou...	11371902	1
Caelo	Dandelion herb withou...	11371904	2
Caelo	Dandelion herb withou...	7031305	1
Caelo	Dandelion herb withou...	11371907	1
Caelo	Dandelion herb withou...	1137907	1
Caelo	Dandelion herb withou...	11217107	1
Caelo	Dandelion herb withou...	13091305	1
Caelo	Dandelion herb withou...	11030311	1
Caelo	Dandelion herb withou...	11217109	1
Caelo	Dandelion root, cut	16161308	1
Caelo	Dandelion root, cut	11003307	1
Caelo	Dandelion root, cut	11003311	2
Caelo	Dandelion root, cut	12363102	2
Caelo	Dandelion root, cut	10024710	1
Caelo	Dandelion root, cut	11003313	1
Caelo	Dandelion root, cut	12363104	3
Caelo	Dandelion root, cut	12363109	2
Caelo	Dandelion root, cut	12363106	1
Caelo	Dandelion root, cut	16161303	1
Caelo	Dandelion root, cut	1140317	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dandelion root, cut	15407703	1
Caelo	Dandelion root, cut	15407703	1
Caelo	Dandelion root, cut	12363103	1
Caelo	Dandelion root, cut	14376602	1
Caelo	Dandelion root, cut	14166504	1
Caelo	Dandelion root, cut	14166502	1
Caelo	Dandelion root, cut	14376603	1
Caelo	Dandelion root, cut	14166501	1
Caelo	Dandelion root, cut	12363107	1
Caelo	Dandelion root, cut	12367916	1
Caelo	Dandelion root, cut	13329702	8
Caelo	Dandelion root, cut	12363108	2
Caelo	Dandelion root, cut	13329705	1
Caelo	Dandelion root, cut	12363110	1
Caelo	Dequalinium chloride	12077510	7
Caelo	Dequalinium chloride	12077508	3
Caelo	Dequalinium chloride	11166203	1
Caelo	Dequalinium chloride	220812K2	1
Caelo	Dequalinium chloride	10879409	1
Caelo	Dequalinium chloride	11166209	1
Caelo	Dequalinium chloride	18079403	1
Caelo	Dequalinium chloride	18079402	1
Caelo	Dequalinium chloride	Ch.B. :311020-6P	1
Caelo	Dequalinium chloride	5315I-01325	1
Caelo	Dequalinium chloride	12361906	2
Caelo	Dequalinium chloride	1405012-02	3
Caelo	Dequalinium chloride	15B16-B04-326872	2
Caelo	Dequalinium chloride	15402801	3
Caelo	Dequalinium chloride	15218406	6
Caelo	Dequalinium chloride	13/26-N10	1
Caelo	Dequalinium chloride	15218404	3
Caelo	Dequalinium chloride	215I-01325	1
Caelo	Dequalinium chloride	15218402	1
Caelo	Dequalinium chloride	14C10-B03	1
Caelo	Dequalinium chloride	18015802	1
Caelo	Dequalinium chloride	14072901	1
Caelo	Dequalinium chloride	12361908	2
Caelo	Dequalinium chloride	18015801	1
Caelo	Dequalinium chloride	15218401	7
Caelo	Dequalinium chloride	15218410	1
Caelo	Dequalinium chloride	13A23-N02	1
Caelo	Dequalinium chloride	15402803	1
Caelo	Dequalinium chloride	15402807	1
Caelo	Dequalinium chloride	12361911	2
Caelo	Dequalinium chloride	12361906	1
Caelo	Dequalinium chloride	15402805	2
Caelo	Dequalinium chloride	3613A-01325	1
Caelo	Dequalinium chloride	152184	1
Caelo	Dequalinium chloride	15218407	1
Caelo	Dequalinium chloride	12361911	1
Caelo	Dequalinium chloride	13348805	2
Caelo	Dequalinium chloride	15B16-B04-342302	2
Caelo	Devil's claw root	6331B110426	2
Caelo	Devil's claw root	12031811	1
Caelo	Devil's claw root	12031806	1
Caelo	Devil's claw root	12031810	1
Caelo	Devil's claw root	12031803	1
Caelo	Devil's claw root	12031809	1
Caelo	Dexamethasone	12165803	1
Caelo	Dexamethasone	13342310	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dexamethasone	15215601	1
Caelo	Dexamethasone	17332414	7
Caelo	Dexamethasone	17332407	2
Caelo	Dexamethasone	15054618	3
Caelo	Dexamethasone	15054617	12
Caelo	Dexamethasone	17332413	2
Caelo	Dexamethasone	1504610	1
Caelo	Dexamethasone	15054614	4
Caelo	Dexamethasone	17332415	1
Caelo	Dexamethasone	17332419	9
Caelo	Dexamethasone	15A05-B01-304637	2
Caelo	Dexamethasone	10369921	2
Caelo	Dexamethasone	17332411	1
Caelo	Dexamethasone	17332415-D20180115Ma	1
Caelo	Dexamethasone	17332416	2
Caelo	Dexamethasone	19022101	5
Caelo	Dexamethasone	17332418	1
Caelo	Dexamethasone	10369920Nachprüfung	1
Caelo	Dexamethasone	1708004-01,180414D	1
Caelo	Dexamethasone	14196402	3
Caelo	Dexamethasone	1335213	1
Caelo	Dexamethasone	15054606	14
Caelo	Dexamethasone	15054605	5
Caelo	Dexamethasone	15054610	2
Caelo	Dexamethasone	15054611	4
Caelo	Dexamethasone	15054609	3
Caelo	Dexamethasone	12033308	3
Caelo	Dexamethasone	6(19086504)	1
Caelo	Dexamethasone	19022107	1
Caelo	Dexamethasone	19022104	1
Caelo	Dexamethasone	19022103	2
Caelo	Dexamethasone	15054603	3
Caelo	Dexamethasone	13352315	2
Caelo	Dexamethasone	160821D	1
Caelo	Dexamethasone	17332412	6
Caelo	Dexamethasone	17104212	1
Caelo	Dexamethasone	17332402	8
Caelo	Dexamethasone	171042165511	1
Caelo	Dexamethasone	171042165585	1
Caelo	Dexamethasone	17104216	4
Caelo	Dexamethasone	1714203	1
Caelo	Dexamethasone	17104210	9
Caelo	Dexamethasone	17104211	16
Caelo	Dexamethasone	17081716	1
Caelo	Dexamethasone	12033301	3
Caelo	Dexamethasone	1704204	1
Caelo	Dexamethasone	17104205	5
Caelo	Dexamethasone	161001-01	1
Caelo	Dexamethasone	1704203	1
Caelo	Dexamethasone	16012207	1
Caelo	Dexamethasone	17104202	1
Caelo	Dexamethasone	1504631	1
Caelo	Dexamethasone	15054616	4
Caelo	Dexamethasone	15054623	5
Caelo	Dexamethasone	9383631	1
Caelo	Dexamethasone	15054621	8
Caelo	Dexamethasone	15054628	10
Caelo	Dexamethasone	15054620/170314D	1
Caelo	Dexamethasone	17332408	6
Caelo	Dexamethasone	17332406	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dexamethasone	15054631	5
Caelo	Dexamethasone	1708804-01	1
Caelo	Dexamethasone	15054627	2
Caelo	Dexamethasone	12033305	4
Caelo	Dexamethasone	1611001-01,170620DM	1
Caelo	Dexamethasone	17104204	7
Caelo	Dexamethasone	10369920	8
Caelo	Dexamethasone	13352304	1
Caelo	Dexamethasone	17332403	14
Caelo	Dexamethasone	13352308	4
Caelo	Dexamethasone	13084305	4
Caelo	Dexamethasone	13084319	5
Caelo	Dexamethasone	13352302	8
Caelo	Dexamethasone	13084320	2
Caelo	Dexamethasone	13084317	6
Caelo	Dexamethasone	13084316	3
Caelo	Dexamethasone	13084309	1
Caelo	Dexamethasone	130854306	1
Caelo	Dexamethasone	13084306	4
Caelo	Dexamethasone	13084307	2
Caelo	Dexamethasone	13084304	2
Caelo	Dexamethasone	13084308	3
Caelo	Dexamethasone	13352312	12
Caelo	Dexamethasone	3084310	1
Caelo	Dexamethasone	13084310	10
Caelo	Dexamethasone	12033317	5
Caelo	Dexamethasone	12033315	5
Caelo	Dexamethasone	12091306	1
Caelo	Dexamethasone	12033318	4
Caelo	Dexamethasone	12033314	8
Caelo	Dexamethasone	12033304	1
Caelo	Dexamethasone	13352310	12
Caelo	Dexamethasone	12033312	2
Caelo	Dexamethasone	13352309	3
Caelo	Dexamethasone acetate	1610662	1
Caelo	Dexamethasone acetate	16106619	1
Caelo	Dexamethasone acetate	16106601	1
Caelo	Dexamethasone acetate	16106602	3
Caelo	Dexamethasone acetate	18173002	5
Caelo	Dexamethasone acetate	16106615	1
Caelo	Dexamethasone acetate	50(18173001)	1
Caelo	Dexamethasone acetate	16106604	1
Caelo	Dexamethasone acetate	14379608	1
Caelo	Dexamethasone acetate	16106608	3
Caelo	Dexamethasone acetate	16106617	2
Caelo	Dexamethasone acetate	16106603	4
Caelo	Dexamethasone acetate	16106611	1
Caelo	Dexamethasone acetate	16106616	8
Caelo	Dexamethasone acetate	16106612	2
Caelo	Dexamethasone acetate	18173012	1
Caelo	Dexamethasone acetate	16106610	1
Caelo	Dexamethasone acetate	19119021	1
Caelo	Dexamethasone acetate	16106603/171025DM	1
Caelo	Dexamethasone acetate	14379615	4
Caelo	Dexamethasone acetate	16106613	5
Caelo	Dexamethasone acetate	16106607	4
Caelo	Dexamethasone acetate	14379601	1
Caelo	Dexamethasone acetate	140502403	1
Caelo	Dexamethasone acetate	14379616	4
Caelo	Dexamethasone acetate	5/AD1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dexamethasone acetate	14200201	2
Caelo	Dexamethasone acetate	13260701	3
Caelo	Dexamethasone acetate	13100403	5
Caelo	Dexamethasone acetate	13260705	1
Caelo	Dexamethasone acetate	3913Q-01326	1
Caelo	Dexamethasone acetate	12021814	2
Caelo	Dexamethasone acetate	13260703	3
Caelo	Dexamethasone acetate	22/AKD2	1
Caelo	Dexamethasone acetate	4113I-01326	1
Caelo	Dexamethasone acetate	13100401	1
Caelo	Dexamethasone acetate	13100404	2
Caelo	Dexamethasone acetate	13100406	2
Caelo	Dexamethasone acetate	13100407	4
Caelo	Dexamethasone acetate	13100408	2
Caelo	Dexamethasone acetate	15243201	1
Caelo	Dexamethasone acetate	12021801	2
Caelo	Dexamethasone acetate	11093215	1
Caelo	Dexamethasone acetate	14379607	1
Caelo	Dexamethasone acetate	14379611	2
Caelo	Dexamethasone acetate	14379609	3
Caelo	Dexamethasone acetate	14379602	2
Caelo	Dexamethasone acetate	15243202	3
Caelo	Dexamethasone acetate	14200208	1
Caelo	Dexamethasone acetate	12021811	1
Caelo	Dexamethasone acetate	14200207	1
Caelo	Dexamethasone acetate	14200205	2
Caelo	Dexamethasone acetate	1405024-01	2
Caelo	Dexamethasone acetate	16106623	1
Caelo	Dexamethasone acetate	12021803	3
Caelo	Dexamethasone acetate	11093211	1
Caelo	Dexamethasone acetate	11093210	1
Caelo	Dexamethasone acetate	13260702	1
Caelo	Dexamethasone acetate	161066165624	1
Caelo	Diclofenac sodium	18039203	1
Caelo	Diclofenac sodium	13214304	1
Caelo	Diclofenac sodium	12214206	1
Caelo	Diclofenac sodium	12214219	1
Caelo	Diclofenac sodium	13214303	1
Caelo	Diclofenac sodium	13214309	3
Caelo	Diclofenac sodium	12214202	1
Caelo	Diclofenac sodium	13214308	2
Caelo	Diclofenac sodium	13214311	9
Caelo	Diclofenac sodium	13214312	2
Caelo	Diclofenac sodium	13214317	5
Caelo	Diclofenac sodium	11091201	1
Caelo	Diclofenac sodium	13214316	1
Caelo	Diclofenac sodium	13214318	3
Caelo	Diclofenac sodium	12L28-B05	1
Caelo	Diclofenac sodium	15169604	6
Caelo	Diclofenac sodium	16349207	1
Caelo	Diclofenac sodium	15278501	1
Caelo	Diclofenac sodium	15278507	1
Caelo	Diclofenac sodium	15278506	6
Caelo	Diclofenac sodium	1503014-01	13
Caelo	Diclofenac sodium	15169613	1
Caelo	Diclofenac sodium	15278505	2
Caelo	Diclofenac sodium	15K23-B03-325753	1
Caelo	Diclofenac sodium	15169675	1
Caelo	Diclofenac sodium	15169615	1
Caelo	Diclofenac sodium	15169609	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Diclofenac sodium	12214205	1
Caelo	Dill fruit	16040406	1
Caelo	Dill fruit	12342908	2
Caelo	Dill fruit	12342903	1
Caelo	Dill fruit	12342906	2
Caelo	Dill fruit	16348702	2
Caelo	Dill fruit	12342909	5
Caelo	Dill fruit	12038302	2
Caelo	Dill fruit	12342904	7
Caelo	Dill fruit	12342901	2
Caelo	Dill fruit	12289901	1
Caelo	Dill fruit	16348707	1
Caelo	Dill fruit	12038304	4
Caelo	Dill fruit	9189003	1
Caelo	Dill herb, cut	13007601	1
Caelo	Dill herb, cut	13007605	1
Caelo	Dill herb, cut	13007602	1
Caelo	Dill herb, cut	11295109	1
Caelo	Dill herb, cut	15178603	1
Caelo	Dill herb, cut	13007608	1
Caelo	Diltiazem hydrochloride	17332021	1
Caelo	Diltiazem hydrochloride	17337009	19
Caelo	Diltiazem hydrochloride	17337016	4
Caelo	Diltiazem hydrochloride	17337018	8
Caelo	Diltiazem hydrochloride	1733712	1
Caelo	Diltiazem hydrochloride	18H09-B03-352908	1
Caelo	Diltiazem hydrochloride	17043213	1
Caelo	Diltiazem hydrochloride	12176804	1
Caelo	Diltiazem hydrochloride	17043212	1
Caelo	Diltiazem hydrochloride	17337012	10
Caelo	Diltiazem hydrochloride	17053181	1
Caelo	Diltiazem hydrochloride	17043207	1
Caelo	Diltiazem hydrochloride	1090917/17043203	1
Caelo	Diltiazem hydrochloride	16190818	1
Caelo	Diltiazem hydrochloride	17043215	11
Caelo	Diltiazem hydrochloride	17043216	15
Caelo	Diltiazem hydrochloride	17043219	17
Caelo	Diltiazem hydrochloride	16190815	1
Caelo	Diltiazem hydrochloride	17043217	18
Caelo	Diltiazem hydrochloride	17043218	9
Caelo	Diltiazem hydrochloride	17337014	3
Caelo	Diltiazem hydrochloride	17043228	6
Caelo	Diltiazem hydrochloride	17h25-b01340699	1
Caelo	Diltiazem hydrochloride	17043224	3
Caelo	Diltiazem hydrochloride	180320DHC/17043218	1
Caelo	Diltiazem hydrochloride	17043223	4
Caelo	Diltiazem hydrochloride	17337007	2
Caelo	Diltiazem hydrochloride	17043230	3
Caelo	Diltiazem hydrochloride	17337004	5
Caelo	Diltiazem hydrochloride	17337010	5
Caelo	Diltiazem hydrochloride	17337011	2
Caelo	Diltiazem hydrochloride	19084106	1
Caelo	Diltiazem hydrochloride	17337025	7
Caelo	Diltiazem hydrochloride	13346905	9
Caelo	Diltiazem hydrochloride	13j24-b04-286571	1
Caelo	Diltiazem hydrochloride	13346904	5
Caelo	Diltiazem hydrochloride	1334609	1
Caelo	Diltiazem hydrochloride	13346901	10
Caelo	Diltiazem hydrochloride	13G16-No5	1
Caelo	Diltiazem hydrochloride	11134113	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Diltiazem hydrochloride	13196103	2
Caelo	Diltiazem hydrochloride	1396103	1
Caelo	Diltiazem hydrochloride	13196105	16
Caelo	Diltiazem hydrochloride	13196107	3
Caelo	Diltiazem hydrochloride	11091301	1
Caelo	Diltiazem hydrochloride	13196104	7
Caelo	Diltiazem hydrochloride	1733701	1
Caelo	Diltiazem hydrochloride	1319604	1
Caelo	Diltiazem hydrochloride	DIL/M-19911	1
Caelo	Diltiazem hydrochloride	3313U-03562	1
Caelo	Diltiazem hydrochloride	10212404	1
Caelo	Diltiazem hydrochloride	12133023	1
Caelo	Diltiazem hydrochloride	12133029	4
Caelo	Diltiazem hydrochloride	170432	1
Caelo	Diltiazem hydrochloride	19084110	1
Caelo	Diltiazem hydrochloride	19084103	1
Caelo	Diltiazem hydrochloride	19084104	2
Caelo	Diltiazem hydrochloride	18301110	3
Caelo	Diltiazem hydrochloride	18301111	6
Caelo	Diltiazem hydrochloride	10212413	1
Caelo	Diltiazem hydrochloride	17337027	7
Caelo	Diltiazem hydrochloride	13346902	6
Caelo	Diltiazem hydrochloride	17043204	8
Caelo	Diltiazem hydrochloride	173327027	1
Caelo	Diltiazem hydrochloride	17043205	5
Caelo	Diltiazem hydrochloride	14059513	17
Caelo	Diltiazem hydrochloride	14059503	3
Caelo	Diltiazem hydrochloride	23121402	1
Caelo	Diltiazem hydrochloride	14059515	2
Caelo	Diltiazem hydrochloride	14318403	4
Caelo	Diltiazem hydrochloride	14318402	4
Caelo	Diltiazem hydrochloride	13031501	1
Caelo	Diltiazem hydrochloride	140595143	1
Caelo	Diltiazem hydrochloride	14318409	14
Caelo	Diltiazem hydrochloride	14059521	2
Caelo	Diltiazem hydrochloride	14C12-B04-294135	2
Caelo	Diltiazem hydrochloride	14318410	3
Caelo	Diltiazem hydrochloride	14318412	19
Caelo	Diltiazem hydrochloride	18061511	1
Caelo	Diltiazem hydrochloride	14318413	3
Caelo	Diltiazem hydrochloride	14318416	2
Caelo	Diltiazem hydrochloride	17043203	12
Caelo	Diltiazem hydrochloride	15/246/DI	1
Caelo	Diltiazem hydrochloride	15212801	8
Caelo	Diltiazem hydrochloride	15212802	7
Caelo	Diltiazem hydrochloride	14059133	1
Caelo	Diltiazem hydrochloride	1091503	1
Caelo	Diltiazem hydrochloride	14059514	18
Caelo	Diltiazem hydrochloride	14059504	3
Caelo	Diltiazem hydrochloride		2
Caelo	Diltiazem hydrochloride	DIL/M-04511	1
Caelo	Diltiazem hydrochloride	13346903	9
Caelo	Diltiazem hydrochloride	25041401	1
Caelo	Diltiazem hydrochloride	12133011	3
Caelo	Diltiazem hydrochloride	6031302	1
Caelo	Diltiazem hydrochloride	12133008	3
Caelo	Diltiazem hydrochloride	12133014	2
Caelo	Diltiazem hydrochloride	12133022	12
Caelo	Diltiazem hydrochloride	1405904	1
Caelo	Diltiazem hydrochloride	10069509	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Diltiazem hydrochloride	2310213562	1
Caelo	Diltiazem hydrochloride	13346908	1
Caelo	Diltiazem hydrochloride	12176801	2
Caelo	Diltiazem hydrochloride	12177101	13
Caelo	Diltiazem hydrochloride	12133006	1
Caelo	Diltiazem hydrochloride	11L20-N06	1
Caelo	Diltiazem hydrochloride	11134117	4
Caelo	Diltiazem hydrochloride	13J24-B04-289432	3
Caelo	Diltiazem hydrochloride	14059502	1
Caelo	Diltiazem hydrochloride	12133028	4
Caelo	Diltiazem hydrochloride	14318414	2
Caelo	Diltiazem hydrochloride	14g22b08299781	1
Caelo	Diltiazem hydrochloride	15230503	11
Caelo	Diltiazem hydrochloride	15414318	1
Caelo	Diltiazem hydrochloride	15414319	23
Caelo	Diltiazem hydrochloride	15414321	2
Caelo	Diltiazem hydrochloride	261115DHC	1
Caelo	Diltiazem hydrochloride	160415D	1
Caelo	Diltiazem hydrochloride	16A18-B07-318733	4
Caelo	Diltiazem hydrochloride	15230513	1
Caelo	Diltiazem hydrochloride	140595	1
Caelo	Diltiazem hydrochloride	1521806	1
Caelo	Diltiazem hydrochloride	15414313	13
Caelo	Diltiazem hydrochloride	15414309	2
Caelo	Diltiazem hydrochloride	15414305	4
Caelo	Diltiazem hydrochloride	15414301	16
Caelo	Diltiazem hydrochloride	16A18-B07-31667	1
Caelo	Diltiazem hydrochloride	15414302	16
Caelo	Diltiazem hydrochloride	15230512	13
Caelo	Diltiazem hydrochloride	15230504	3
Caelo	Diltiazem hydrochloride	16A18-B07-319570	2
Caelo	Diltiazem hydrochloride	15230516	1
Caelo	Diltiazem hydrochloride	15414303	1
Caelo	Diltiazem hydrochloride	154114319	1
Caelo	Diltiazem hydrochloride	170729181	1
Caelo	Diltiazem hydrochloride	17043202	16
Caelo	Diltiazem hydrochloride	1020617/16190813	1
Caelo	Diltiazem hydrochloride	170602DH/16190813	1
Caelo	Diltiazem hydrochloride	16190814	7
Caelo	Diltiazem hydrochloride	16190817	1
Caelo	Diltiazem hydrochloride	1250417	1
Caelo	Diltiazem hydrochloride	15435807	1
Caelo	Diltiazem hydrochloride	15230502	12
Caelo	Diltiazem hydrochloride	16190801	1
Caelo	Diltiazem hydrochloride	16190813	14
Caelo	Diltiazem hydrochloride	15414312	1
Caelo	Diltiazem hydrochloride	16190805	3
Caelo	Diltiazem hydrochloride	16190807	3
Caelo	Diltiazem hydrochloride	15414310	1
Caelo	Diltiazem hydrochloride	16190810	3
Caelo	Diltiazem hydrochloride	15414320	1
Caelo	Diltiazem hydrochloride	16190806	16
Caelo	Diltiazem hydrochloride	161908	1
Caelo	Diltiazem hydrochloride	15230507	5
Caelo	Diltiazem hydrochloride	15230515	1
Caelo	Diltiazem hydrochloride	15212806	2
Caelo	Diltiazem hydrochloride	15230501	3
Caelo	Diltiazem hydrochloride	15230000000000000000	1
Caelo	Diltiazem hydrochloride	15230510	3
Caelo	Diphenhydramine hydro...	12007716	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Diphenhydramine hydro...	2215I-01339	1
Caelo	Diphenhydramine hydro...	15220404	3
Caelo	Diphenhydramine hydro...	15220403	3
Caelo	Diphenhydramine hydro...	15220402	2
Caelo	Diphenhydramine hydro...	15072603	5
Caelo	Diphenhydramine hydro...	15072630	1
Caelo	Diphenhydramine hydro...	15072602	3
Caelo	Diphenhydramine hydro...	13331813	1
Caelo	Diphenhydramine hydro...	20071502	1
Caelo	Diphenhydramine hydro...	13331810	1
Caelo	Diphenhydramine hydro...	1814A-01339	1
Caelo	Diphenhydramine hydro...	12329815	2
Caelo	Diphenhydramine hydro...	13343104	1
Caelo	Diphenhydramine hydro...	15220405	2
Caelo	Diphenhydramine hydro...	1912I-01339	1
Caelo	Diphenhydramine hydro...	2041403	1
Caelo	Diphenhydramine hydro...	12007721	2
Caelo	Diphenhydramine hydro...	12007715	2
Caelo	Diphenhydramine hydro...	12007713	2
Caelo	Diphenhydramine hydro...	25071306	1
Caelo	Diphenhydramine hydro...	12007702	1
Caelo	Diphenhydramine hydro...	5121204	1
Caelo	Diphenhydramine hydro...	10122030	1
Caelo	Diphenhydramine hydro...	31031508	1
Caelo	Diphenhydramine hydro...	15072608	1
Caelo	Diphenhydramine hydro...	13100462	1
Caelo	Diphenhydramine hydro...	15072614	7
Caelo	Diphenhydramine hydro...	3071715	1
Caelo	Diphenhydramine hydro...	16004908	2
Caelo	Diphenhydramine hydro...	16004910	4
Caelo	Diphenhydramine hydro...	18008201	3
Caelo	Diphenhydramine hydro...	18008205	1
Caelo	Diphenhydramine hydro...	15072609	4
Caelo	Dithranol	13120504	1
Caelo	Dithranol	1212I-01274	1
Caelo	Dithranol	12375305	1
Caelo	Dithranol	12375601	1
Caelo	Dithranol	12375403	1
Caelo	Dithranol	12375302	4
Caelo	Dithranol	120145041.2013	1
Caelo	Dithranol	12014502	1
Caelo	Dithranol	27111206	1
Caelo	Dithranol	13120503	1
Caelo	Dithranol	APDT/3652011D	1
Caelo	Dithranol	12014504	2
Caelo	Dithranol	11064501	1
Caelo	Dithranol	12236811	4
Caelo	Dithranol	18115306	1
Caelo	Dithranol	12236809	1
Caelo	Dithranol	18B23-B05-348137	1
Caelo	Dithranol	13120502	1
Caelo	Dithranol	18115302	1
Caelo	Dithranol	13120702	1
Caelo	Dithranol	14101715	1
Caelo	Doxycycline hyclate	16I06-B02-341857	1
Caelo	Doxycycline hyclate	16I06-B02-335675	2
Caelo	Doxylamine hydrogen s...	14k27-B07	1
Caelo	Doxylamine hydrogen s...	15F01-B04-309189	3
Caelo	Doxylamine hydrogen s...	15K04B01327356	1
Caelo	Doxylamine hydrogen s...	15F01-B04-317385	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Dry beads Caelo ora...	11005811	1
Caelo	Echinacea pallida roo...	12244504	1
Caelo	Echinacea pallida roo...	11115206	1
Caelo	Econazole nitrate	ECN140301	1
Caelo	Elder berries	12267801	1
Caelo	Elder flower, rubbed	13244601	4
Caelo	Elder flower, rubbed	13244602	2
Caelo	Elder flower, rubbed	12221904	10
Caelo	Elder flower, rubbed	31051301	1
Caelo	Elder flower, rubbed	12221909	1
Caelo	Elder flower, rubbed	5122013d	1
Caelo	Elder flower, rubbed	14111317	1
Caelo	Elder flower, rubbed	12221910	2
Caelo	Elder flower, rubbed	12221906	1
Caelo	Elder flower, rubbed	12221907	4
Caelo	Elder flower, rubbed	12221905	7
Caelo	Elder flower, rubbed	10291310	5
Caelo	Elder flower, rubbed	1029310	1
Caelo	Elder flower, rubbed	1571A121022	14
Caelo	Elder flower, rubbed	12221904-3360319071	1
Caelo	Elder flower, rubbed	10291312	2
Caelo	Elder flower, rubbed	10291308	2
Caelo	Elder flower, rubbed	10291309	3
Caelo	Elder flower, rubbed	12221902	10
Caelo	Elder flower, rubbed	1571a121022	1
Caelo	Elecampane rhizome, cut	12113304	3
Caelo	Elecampane rhizome, cut	1211334	1
Caelo	Elecampane rhizome, cut	13152401	1
Caelo	Elecampane rhizome, cut	12113303	1
Caelo	Elecampane rhizome, cut	13152402	1
Caelo	Elecampane rhizome, cut	13152403	2
Caelo	Elecampane rhizome, cut	11374703	1
Caelo	Elecampane rhizome, cut	11374702	1
Caelo	Elecampane rhizome, cut	6341B	1
Caelo	Eleutherococcus, cut	10381817	1
Caelo	Eleutherococcus, cut	12326701	1
Caelo	Eleutherococcus, cut	12326707	1
Caelo	Eleutherococcus, cut	12326703	1
Caelo	Eleutherococcus, cut	10381818	2
Caelo	Elm bark, cut	10347204	1
Caelo	Elm bark, cut	12287902	2
Caelo	Elm bark, cut	10347203	1
Caelo	Emulsifying cetostear...	16238607	9
Caelo	Emulsifying cetostear...	18018006	1
Caelo	Emulsifying cetostear...	15385602	1
Caelo	Emulsifying cetostear...	31032015B	1
Caelo	Emulsifying cetostear...	15385602	2
Caelo	Emulsifying cetostear...	13233006	1
Caelo	Emulsifying cetostear...	16119607	2
Caelo	Emulsifying cetostear...	17056802	1
Caelo	Emulsifying cetostear...	13144105	1
Caelo	Emulsifying cetostear...	16I13-B01-327864	1
Caelo	Emulsifying cetostear...	608/10072017	1
Caelo	Emulsifying cetostear...	13352102	1
Caelo	Emulsifying cetostear...	15071704	1
Caelo	Emulsifying cetostear...	14322502	2
Caelo	Emulsifying cetostear...	5011505	1
Caelo	Emulsifying cetostear...	14158202	1
Caelo	Emulsifying cetostear...	14322505	7
Caelo	Emulsifying cetostear...	14322505	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Emulsifying cetostear...	15271408	3
Caelo	Emulsifying cetostear...	23051503	1
Caelo	Emulsifying cetostear...	15271402	9
Caelo	Emulsifying cetostear...	15271403	10
Caelo	Emulsifying cetostear...	15271405	10
Caelo	Emulsifying cetostear...	15090706	3
Caelo	Emulsifying cetostear...	15090703	2
Caelo	Emulsifying cetostear...	14322512	1
Caelo	Emulsifying cetostear...	14156201	1
Caelo	Emulsifying cetostear...	14G31-B03-305612	3
Caelo	Emulsifying cetostear...	14322506	1
Caelo	Emulsifying cetostear...	27/AKD\$	1
Caelo	Emulsifying cetostear...	15090701	9
Caelo	Emulsifying cetostear...	12331402	1
Caelo	Emulsifying cetostear...	15271410	1
Caelo	Emulsifying cetostear...	15385606	1
Caelo	Emulsifying cetostear...	13144103	1
Caelo	Emulsifying cetostear...	17056806	5
Caelo	Emulsifying cetostear...	13233002	1
Caelo	Emulsifying cetostear...	9170717	1
Caelo	Emulsifying cetostear...	14158201	6
Caelo	Emulsifying cetostear...	13144107	2
Caelo	Emulsifying cetostear...	12134203	1
Caelo	Emulsifying cetostear...	12331407	3
Caelo	Emulsifying cetostear...	25/AKD1	1
Caelo	Emulsifying cetostear...	13144104	3
Caelo	Emulsifying cetostear...	12331402	5
Caelo	Emulsifying cetostear...	12231402	1
Caelo	Emulsifying cetostear...	11291103	1
Caelo	Emulsifying cetostear...	12229402	2
Caelo	Emulsifying cetostear...	12065709	6
Caelo	Emulsifying cetostear...	11009204	1
Caelo	Emulsifying cetostear...	12C30-N01	2
Caelo	Emulsifying cetostear...	18081201	1
Caelo	Emulsifying cetostear...	25051201	1
Caelo	Emulsifying cetostear...	11291101	1
Caelo	Emulsifying cetostear...	17139201	4
Caelo	Emulsifying cetostear...	13144105	11
Caelo	Emulsifying cetostear...	24101305	1
Caelo	Emulsifying cetostear...	2011404	1
Caelo	Emulsifying cetostear...	17198003	16
Caelo	Emulsifying cetostear...	13233005	1
Caelo	Emulsifying cetostear...	13233008	3
Caelo	Emulsifying cetostear...	16238604	1
Caelo	Emulsifying cetostear...	14158201	4
Caelo	Emulsifying cetostear...	17261701	1
Caelo	Emulsifying cetostear...	13233001	1
Caelo	Emulsifying cetostear...	13L17-B05-289806	1
Caelo	Emulsifying cetostear...	23091402	1
Caelo	Emulsifying cetostear...	17261702	2
Caelo	Emulsifying cetostear...	Unbekannt	1
Caelo	Emulsifying cetostear...	17198002	3
Caelo	Emulsifying cetostear...	29091803	1
Caelo	Emulsifying cetostear...	18018005	2
Caelo	Emulsifying cetostear...	19032014B	1
Caelo	Emulsifying cetostear...	13233002	1
Caelo	Emulsifying cetostear...	15271407	6
Caelo	Eosin Y	12e25n01	1
Caelo	Eosin Y	13I11-N03	1
Caelo	Eosin Y	14H28-M05	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Eosin Y	13F14-T01-017395	7
Caelo	Eosin Y	12e25-n01	1
Caelo	Ephedrine hydrochloride	11126407	1
Caelo	Ephedrine hydrochloride	10103425	1
Caelo	Ephedrine hydrochloride	16119715	2
Caelo	Ephedrine hydrochloride	18006701	1
Caelo	Ephedrine hydrochloride	16119715200218P301	1
Caelo	Ephedrine hydrochloride	16119709	1
Caelo	Ephedrine hydrochloride	10119709	1
Caelo	Ephedrine hydrochloride	1619704	1
Caelo	Ephedrine hydrochloride	16119702	1
Caelo	Ephedrine hydrochloride	15179302	1
Caelo	Ephedrine hydrochloride	16119703	1
Caelo	Ephedrine hydrochloride	14378507	1
Caelo	Ephedrine hydrochloride	151J9304	1
Caelo	Ephedrine hydrochloride	18006709	1
Caelo	Ephedrine hydrochloride	14378506	2
Caelo	Ephedrine hydrochloride	14378502	1
Caelo	Ephedrine hydrochloride	13323316	1
Caelo	Ephedrine hydrochloride	13323310	1
Caelo	Ephedrine hydrochloride	13323305	2
Caelo	Ephedrine hydrochloride	13323302	2
Caelo	Ephedrine hydrochloride	12017223	1
Caelo	Ephedrine hydrochloride	12017217	1
Caelo	Ephedrine hydrochloride	12017215	1
Caelo	Ephedrine hydrochloride	14378510	1
Caelo	Ephedrine hydrochloride	12017209	2
Caelo	Epinephrine hydrogen ...	14301002	2
Caelo	Epinephrine hydrogen ...	14007808	1
Caelo	Epinephrine hydrogen ...	12337801	1
Caelo	Epinephrine hydrogen ...	12120202	1
Caelo	Epinephrine hydrogen ...	12337805	3
Caelo	Epinephrine hydrogen ...	10237109	1
Caelo	Epinephrine hydrogen ...	15194102	1
Caelo	Epinephrine hydrogen ...	123337805	1
Caelo	Epinephrine hydrogen ...	15239703	1
Caelo	Epinephrine hydrogen ...	15239701	2
Caelo	Epinephrine hydrogen ...	210316-3	1
Caelo	Epinephrine hydrogen ...	15194104	1
Caelo	Epinephrine hydrogen ...	16276402	1
Caelo	Epinephrine hydrogen ...	162746402	1
Caelo	Epinephrine hydrogen ...	17111503	2
Caelo	Epinephrine hydrogen ...	16276406	1
Caelo	Epinephrine hydrogen ...	18179503	1
Caelo	Epinephrine hydrogen ...	18179504	1
Caelo	Epinephrine hydrogen ...	210316-1	1
Caelo	Equisetum stem, cut	11076307	1
Caelo	Equisetum stem, cut	11076303	13
Caelo	Equisetum stem, cut	20111212	1
Caelo	Equisetum stem, cut	11021904	1
Caelo	Equisetum stem, cut	17101203	1
Caelo	Equisetum stem, cut	8101201	1
Caelo	Equisetum stem, cut		2
Caelo	Equisetum stem, cut	3091201	1
Caelo	Equisetum stem, cut	11076302	2
Caelo	Equisetum stem, cut	11021901	1
Caelo	Equisetum stem, cut	10375507	2
Caelo	Equisetum stem, cut	10328915	1
Caelo	Equisetum stem, cut	12288804	11
Caelo	Equisetum stem, cut	11076306	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Equisetum stem, cut	10375520	1
Caelo	Equisetum stem, cut	13159702	10
Caelo	Equisetum stem, cut	10328901	1
Caelo	Equisetum stem, cut	13350309	1
Caelo	Equisetum stem, cut	13159705-3387310190	1
Caelo	Equisetum stem, cut	30101305	1
Caelo	Equisetum stem, cut	13200201	1
Caelo	Equisetum stem, cut	2658181	1
Caelo	Equisetum stem, cut	13350301	1
Caelo	Equisetum stem, cut	13344302	3
Caelo	Equisetum stem, cut	13159702-3464312101	1
Caelo	Equisetum stem, cut	282573	1
Caelo	Equisetum stem, cut	2011408	1
Caelo	Equisetum stem, cut	13344311	1
Caelo	Equisetum stem, cut	13344306	1
Caelo	Equisetum stem, cut	13344303	6
Caelo	Equisetum stem, cut	13344309	3
Caelo	Equisetum stem, cut	13344305	1
Caelo	Equisetum stem, cut	15015404	1
Caelo	Equisetum stem, cut	9101313	1
Caelo	Equisetum stem, cut	13344304	5
Caelo	Equisetum stem, cut	311042	1
Caelo	Equisetum stem, cut	13153702	1
Caelo	Equisetum stem, cut	13159706	5
Caelo	Equisetum stem, cut	6081304	1
Caelo	Equisetum stem, cut	13091301	1
Caelo	Equisetum stem, cut	12288802	2
Caelo	Equisetum stem, cut	13159705	6
Caelo	Equisetum stem, cut	1359702	1
Caelo	Equisetum stem, cut	12288809	3
Caelo	Equisetum stem, cut	12288808	8
Caelo	Equisetum stem, cut	8051302	1
Caelo	Equisetum stem, cut	8041301	1
Caelo	Equisetum stem, cut	1288804	1
Caelo	Equisetum stem, cut	13159704	6
Caelo	Erythromycin	13175802	8
Caelo	Erythromycin	13175604	7
Caelo	Erythromycin	24101312	1
Caelo	Erythromycin	13145504	1
Caelo	Erythromycin	13175504-3413310132	1
Caelo	Erythromycin	13175804	2
Caelo	Erythromycin	13175506-3559412014	1
Caelo	Erythromycin	13175612	2
Caelo	Erythromycin	59/AKD59	1
Caelo	Erythromycin	13175607	7
Caelo	Erythromycin	132152108	1
Caelo	Erythromycin	13175601	10
Caelo	Erythromycin	20122013D	1
Caelo	Erythromycin	1308016-01	31
Caelo	Erythromycin	13175504	6
Caelo	Erythromycin	13175505	3
Caelo	Erythromycin	131756017	1
Caelo	Erythromycin	511111	1
Caelo	Erythromycin	13175808	2
Caelo	Erythromycin	13175608	11
Caelo	Erythromycin	14044008	4
Caelo	Erythromycin	14044105	4
Caelo	Erythromycin	1404409	1
Caelo	Erythromycin	14044010	4
Caelo	Erythromycin	14044007	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Erythromycin	14044002	7
Caelo	Erythromycin	14b06-b35-293508	1
Caelo	Erythromycin	1414M-01362	1
Caelo	Erythromycin	14044009	2
Caelo	Erythromycin	13175801	9
Caelo	Erythromycin	14044201	5
Caelo	Erythromycin	14044004	12
Caelo	Erythromycin	14044005	7
Caelo	Erythromycin	13175611	5
Caelo	Erythromycin	13175809	5
Caelo	Erythromycin	14044001	10
Caelo	Erythromycin	1114Q-01362	1
Caelo	Erythromycin	13175811	8
Caelo	Erythromycin	13178756tzgklhj	1
Caelo	Erythromycin	13175803	7
Caelo	Erythromycin	3092014B	1
Caelo	Erythromycin	12156203	2
Caelo	Erythromycin	18096201	3
Caelo	Erythromycin	12152107	1
Caelo	Erythromycin	12152101	3
Caelo	Erythromycin	17088309	3
Caelo	Erythromycin	3110195	1
Caelo	Erythromycin	12152205	10
Caelo	Erythromycin	12152202	4
Caelo	Erythromycin	20270912-3	1
Caelo	Erythromycin	1215219	1
Caelo	Erythromycin	12152109	10
Caelo	Erythromycin	22041302	1
Caelo	Erythromycin	11195405	1
Caelo	Erythromycin	11195402	2
Caelo	Erythromycin	1195407	1
Caelo	Erythromycin	12152102	3
Caelo	Erythromycin	11195407	4
Caelo	Erythromycin	11195408	1
Caelo	Erythromycin	27071202	1
Caelo	Erythromycin	11j06-n05	1
Caelo	Erythromycin	14044202	3
Caelo	Erythromycin	12152105	4
Caelo	Erythromycin	13175503	9
Caelo	Erythromycin	12156201	3
Caelo	Erythromycin	12152208	4
Caelo	Erythromycin	13175506	5
Caelo	Erythromycin	13189006	3
Caelo	Erythromycin	13175603	12
Caelo	Erythromycin	13175502	5
Caelo	Erythromycin	13189007	2
Caelo	Erythromycin	13189004	5
Caelo	Erythromycin	1315501	1
Caelo	Erythromycin	13175501	5
Caelo	Erythromycin	12156202	7
Caelo	Erythromycin	13189002	7
Caelo	Erythromycin	12156204	4
Caelo	Erythromycin	13189001	7
Caelo	Erythromycin	11195306	1
Caelo	Erythromycin	12156206	2
Caelo	Erythromycin	19071306	1
Caelo	Erythromycin	12156207	4
Caelo	Erythromycin	12156208	7
Caelo	Erythromycin	12152104	2
Caelo	Erythromycin	1008018-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Erythromycin	14044107	3
Caelo	Erythromycin	16079608	17
Caelo	Erythromycin	16079213	13
Caelo	Erythromycin	17088008-20180206Ma2	1
Caelo	Erythromycin	16079202	6
Caelo	Erythromycin	140930E	1
Caelo	Erythromycin	16079201	6
Caelo	Erythromycin	15042902	3
Caelo	Erythromycin	160703801	8
Caelo	Erythromycin	1504812	1
Caelo	Erythromycin	15042813, 20161112M	1
Caelo	Erythromycin	15042813	12
Caelo	Erythromycin	1504813	1
Caelo	Erythromycin	16079208	8
Caelo	Erythromycin	16079203	17
Caelo	Erythromycin	1181116	1
Caelo	Erythromycin	1171116	1
Caelo	Erythromycin	15042816	4
Caelo	Erythromycin	16979201	1
Caelo	Erythromycin	15042814	10
Caelo	Erythromycin	16079204	11
Caelo	Erythromycin	15042901	7
Caelo	Erythromycin	150428	1
Caelo	Erythromycin	17088013	9
Caelo	Erythromycin	17088001	2
Caelo	Erythromycin	15042809	3
Caelo	Erythromycin	15245005	1
Caelo	Erythromycin	1524004	1
Caelo	Erythromycin	15042804	1
Caelo	Erythromycin	15042806	11
Caelo	Erythromycin	17088406	5
Caelo	Erythromycin	170718A	1
Caelo	Erythromycin	17088403-20180619Ma1	1
Caelo	Erythromycin	178824	1
Caelo	Erythromycin	17088203	9
Caelo	Erythromycin	17088204	20
Caelo	Erythromycin	16079509	5
Caelo	Erythromycin	1788402	1
Caelo	Erythromycin	17088011	1
Caelo	Erythromycin	17088405	7
Caelo	Erythromycin	17088010-20180404Ma2	1
Caelo	Erythromycin	17088403	9
Caelo	Erythromycin	17088401	1
Caelo	Erythromycin	17088301	25
Caelo	Erythromycin	15042812	1
Caelo	Erythromycin	16079210	8
Caelo	Erythromycin	16G06-803	1
Caelo	Erythromycin	16079501	13
Caelo	Erythromycin	17088003	6
Caelo	Erythromycin	16079503	4
Caelo	Erythromycin	17088002	12
Caelo	Erythromycin	17088007	8
Caelo	Erythromycin	16079214-20171127Ma6	1
Caelo	Erythromycin	16079510	4
Caelo	Erythromycin	1708808	1
Caelo	Erythromycin	16079603	6
Caelo	Erythromycin	16079702	6
Caelo	Erythromycin	16079704	2
Caelo	Erythromycin	16079614	5
Caelo	Erythromycin	16079609	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Erythromycin	160796	1
Caelo	Erythromycin	16079612	11
Caelo	Erythromycin	1708807	1
Caelo	Erythromycin	16079607	1
Caelo	Erythromycin	16079605	2
Caelo	Erythromycin	16079604	4
Caelo	Erythromycin	16079514	1
Caelo	Erythromycin	16079214-20170704Ma1	1
Caelo	Erythromycin	15042608	3
Caelo	Erythromycin	16079610	1
Caelo	Erythromycin	17088010	9
Caelo	Erythromycin	16079505	3
Caelo	Erythromycin	16079214	8
Caelo	Erythromycin	16079216	7
Caelo	Erythromycin	16079207	5
Caelo	Erythromycin	16079511	6
Caelo	Erythromycin	15042815	3
Caelo	Erythromycin	150242813	1
Caelo	Erythromycin	16079602	5
Caelo	Erythromycin	16079506	3
Caelo	Erythromycin	15042816, 20170427Ma1	1
Caelo	Erythromycin	16079513	8
Caelo	Erythromycin	16079502	11
Caelo	Erythromycin	17088009	5
Caelo	Erythromycin	16079507	1
Caelo	Erythromycin	16079515	1
Caelo	Erythromycin	17088402	9
Caelo	Erythromycin	16079512	4
Caelo	Erythromycin	17088008	9
Caelo	Erythromycin	17088008-E20180112Ma	1
Caelo	Erythromycin	16G06B03329150	1
Caelo	Erythromycin	15245007	2
Caelo	Erythromycin	1705037-01	1
Caelo	Erythromycin	15042808	1
Caelo	Erythromycin	14/30-B01-300059	1
Caelo	Erythromycin	15042601	2
Caelo	Erythromycin	14044014	2
Caelo	Erythromycin	17088207	1
Caelo	Erythromycin	14044207	1
Caelo	Erythromycin	1404418	1
Caelo	Erythromycin	17088306	13
Caelo	Erythromycin	17088308	2
Caelo	Erythromycin	17088304	1
Caelo	Erythromycin	17088305	5
Caelo	Erythromycin	14044011	2
Caelo	Erythromycin	18096401	1
Caelo	Erythromycin	15042612	6
Caelo	Erythromycin	15245008	4
Caelo	Erythromycin	14/30-B01-308189	4
Caelo	Erythromycin	15042702	4
Caelo	Erythromycin	15042610	12
Caelo	Erythromycin	15042604	3
Caelo	Erythromycin	17088302	2
Caelo	Erythromycin	17088210-E20181115Ma	1
Caelo	Erythromycin	15042602	5
Caelo	Erythromycin	EB06/14	1
Caelo	Erythromycin	18096209	1
Caelo	Erythromycin	18096210	1
Caelo	Erythromycin	14F13-B02-296464	5
Caelo	Erythromycin	14044006	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Erythromycin	14044101	7
Caelo	Erythromycin	14044115	6
Caelo	Erythromycin	14044109	5
Caelo	Erythromycin	140440044	1
Caelo	Erythromycin	5052015D	1
Caelo	Erythromycin	1430b01300059	1
Caelo	Erythromycin	17068309	1
Caelo	Erythromycin	14044206	2
Caelo	Erythromycin	14F13-B02-300061	3
Caelo	Erythromycin	18096202	3
Caelo	Erythromycin	18096208	1
Caelo	Erythromycin	14044114	6
Caelo	Erythromycin	14044106	1
Caelo	Erythromycin	17088404	9
Caelo	Erythromycin	14044118	3
Caelo	Erythromycin	14044012	2
Caelo	Erythromycin	1708/8308	1
Caelo	Erythromycin	18096211	3
Caelo	Erythromycin	180962011905013-04	1
Caelo	Erythromycin	17088403-E2081115Ma1	1
Caelo	Erythromycin	14 30-B01-306927	1
Caelo	Erythromycin	15245101	10
Caelo	Erythromycin	15042705	12
Caelo	Erythromycin	16i24-B10	1
Caelo	Erythromycin	15042609	9
Caelo	Erythromycin	15042704	7
Caelo	Erythromycin	15042606	1
Caelo	Erythromycin	15042701	9
Caelo	Erythromycin	15042611	6
Caelo	Erythromycin	1780210	1
Caelo	Erythromycin	15245004	12
Caelo	Erythromycin	16J24-B-341876	1
Caelo	Erythromycin	1709038-01	1
Caelo	Erythromycin	15042611-20160504Ma1	1
Caelo	Erythromycin	1524002	1
Caelo	Erythromycin	15042711	12
Caelo	Erythromycin	15245003	8
Caelo	Erythromycin	15042803	10
Caelo	Erythromycin	15042802	3
Caelo	Erythromycin	170880103	1
Caelo	Erythromycin	3081601	1
Caelo	Erythromycin	15042706	8
Caelo	Erythromycin	15042707	1
Caelo	Erythromycin	16073038-02	1
Caelo	Erythromycin	17088212	1
Caelo	Erythromycin	15042703	7
Caelo	Erythromycin	1404414	1
Caelo	Erythromycin	15042613	3
Caelo	Estradiol benzoate	17357605	1
Caelo	Estradiol benzoate	18144405	1
Caelo	Estradiol benzoate	15078603	16
Caelo	Estradiol benzoate	6091305	1
Caelo	Estradiol benzoate	15078606	5
Caelo	Estradiol benzoate	10124912	1
Caelo	Estradiol benzoate	12255905	10
Caelo	Estradiol benzoate	10124919	1
Caelo	Estradiol benzoate		3
Caelo	Estradiol benzoate	10124909	4
Caelo	Estradiol benzoate	12255904	10
Caelo	Estradiol benzoate	10124921	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Estradiol benzoate	10124918	2
Caelo	Estradiol benzoate	10124915	2
Caelo	Estradiol benzoate	15078603	1
Caelo	Estradiol benzoate	10124909	1
Caelo	Estradiol benzoate	84303830	1
Caelo	Estradiol benzoate	15C10-B04-305523	1
Caelo	Estradiol benzoate	15078604	16
Caelo	Estradiol benzoate	15078604	3
Caelo	Estradiol benzoate	15C10B04305523	1
Caelo	Estradiol benzoate	1507806	1
Caelo	Estradiol benzoate	92137389	1
Caelo	Estradiol benzoate	15078602	6
Caelo	Estradiol benzoate	15078608	12
Caelo	Estradiol benzoate	13277605	3
Caelo	Estradiol benzoate	16003702	3
Caelo	Estradiol benzoate	13277602	8
Caelo	Estradiol benzoate	ch.b.150700901	1
Caelo	Estradiol benzoate	16003701	4
Caelo	Estradiol benzoate	21121204	1
Caelo	Estradiol benzoate	150918EB	1
Caelo	Estradiol benzoate	1100417	1
Caelo	Estradiol benzoate	15078609	1
Caelo	Estradiol benzoate	1149	1
Caelo	Estradiol benzoate	17357602	4
Caelo	Estradiol benzoate	13276401	6
Caelo	Estradiol benzoate	81403700	1
Caelo	Estradiol benzoate	111004-02	1
Caelo	Estradiol benzoate	18144402	8
Caelo	Estradiol benzoate	13277607	3
Caelo	Estradiol hemihydrate	15385410	1
Caelo	Estradiol hemihydrate	11078105	1
Caelo	Estradiol hemihydrate	15D01-B02-322166	1
Caelo	Estradiol hemihydrate	9248916	2
Caelo	Estradiol hemihydrate	15038014	2
Caelo	Estradiol hemihydrate	84201560	2
Caelo	Estradiol hemihydrate	13285804	1
Caelo	Estradiol hemihydrate	15038010	3
Caelo	Estradiol hemihydrate	1503801	1
Caelo	Estradiol hemihydrate	18082105	3
Caelo	Estradiol hemihydrate	1705018-01	21
Caelo	Estradiol hemihydrate	18082101	1
Caelo	Estradiol hemihydrate	18082104	2
Caelo	Estradiol hemihydrate	12268012	4
Caelo	Estradiol hemihydrate	22101301	1
Caelo	Estradiol hemihydrate	11078113	1
Caelo	Estradiol hemihydrate	11078104	1
Caelo	Estradiol hemihydrate	12172003	2
Caelo	Estradiol hemihydrate	13285802	5
Caelo	Estradiol hemihydrate	13285801	1
Caelo	Estradiol hemihydrate	15385415	2
Caelo	Estradiol hemihydrate	12268003	1
Caelo	Estradiol hemihydrate	13285811	3
Caelo	Estradiol hemihydrate	19051501	1
Caelo	Estradiol valerate	18067403	1
Caelo	Estradiol valerate	16327709	3
Caelo	Estradiol valerate	18067407	1
Caelo	Estradiol valerate	18067409	4
Caelo	Estradiol valerate	16327706	2
Caelo	Estradiol valerate	1810013	1
Caelo	Estradiol valerate	16327703	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Estradiol valerate	16327702	3
Caelo	Estradiol valerate	15c10-b05-314840	1
Caelo	Estradiol valerate	5315A-03259	1
Caelo	Estradiol valerate	16327710	1
Caelo	Estradiol valerate	15c10-B05-305396	1
Caelo	Estradiol valerate	15194808	1
Caelo	Estradiol valerate	15C10B05308851	1
Caelo	Estradiol valerate	18067402	8
Caelo	Estradiol valerate	4214U-03259	1
Caelo	Estradiol valerate	16327711	1
Caelo	Estradiol valerate	14E26-B02-294844	1
Caelo	Estradiol valerate	13264001	1
Caelo	Estradiol valerate	12G05-N09	4
Caelo	Estradiol valerate	18067405	1
Caelo	Estradiol valerate	test	1
Caelo	Estradiol valerate	13264003	1
Caelo	Estradiol valerate	18062402	1
Caelo	Estriol	16053405	4
Caelo	Estriol	15297805	1
Caelo	Estriol	15028905	2
Caelo	Estriol	15008902	1
Caelo	Estriol	12H06N02	1
Caelo	Estriol	14077904	1
Caelo	Estriol	15028904	1
Caelo	Estriol	14j20-b07	1
Caelo	Estriol	16053408	5
Caelo	Estriol	17133704	11
Caelo	Estriol	16138204	4
Caelo	Estriol	19035607	1
Caelo	Estriol	18247101	2
Caelo	Estriol	17348707	2
Caelo	Estriol	17348706	3
Caelo	Estriol	17348702	7
Caelo	Estriol	17348701	7
Caelo	Estriol	17133710	2
Caelo	Estriol	17133709	8
Caelo	Estriol	1713373	1
Caelo	Estriol	1733706	1
Caelo	Estriol	17133703	2
Caelo	Estriol	17133706	17
Caelo	Estriol	16/19-B01-331467	1
Caelo	Estriol	17133705	2
Caelo	Estriol	16053409	1
Caelo	Estriol	16317702	6
Caelo	Estriol	16317703	4
Caelo	Estriol	16053404	2
Caelo	Estriol	12032503	4
Caelo	Estriol	15297804	14
Caelo	Estriol	14272101	2
Caelo	Estriol	27061303	1
Caelo	Estriol	12298707	2
Caelo	Estriol	12298706	1
Caelo	Estriol	12298702	2
Caelo	Estriol	12298704	6
Caelo	Estriol	11254412	1
Caelo	Estriol	13256004	11
Caelo	Estriol	17042014B	1
Caelo	Estriol	11254410	1
Caelo	Estriol	13256001	1
Caelo	Estriol	11254402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Estriol	14077903	2
Caelo	Estriol	14077901	1
Caelo	Estriol	14272101	1
Caelo	Estriol	1409018-01	1
Caelo	Estriol	10053206	1
Caelo	Estriol	15028901	5
Caelo	Estriol	10H30-N06	1
Caelo	Estriol	12032507	5
Caelo	Estriol	14272104	3
Caelo	Estriol	1477904	1
Caelo	Estriol	14077904	10
Caelo	Estriol	24031501	1
Caelo	Ethacridine lactate m...	16273403	7
Caelo	Ethacridine lactate m...	15438009	2
Caelo	Ethacridine lactate m...	15438007	1
Caelo	Ethacridine lactate m...	11J07-N05	1
Caelo	Ethacridine lactate m...	15438003	4
Caelo	Ethacridine lactate m...	141203401	1
Caelo	Ethacridine lactate m...	15177608	1
Caelo	Ethacridine lactate m...	12149507	1
Caelo	Ethacridine lactate m...	15177603	4
Caelo	Ethacridine lactate m...	15177612	1
Caelo	Ethacridine lactate m...	12149507	2
Caelo	Ethacridine lactate m...	10288702	1
Caelo	Ethacridine lactate m...	15438007	4
Caelo	Ethacridine lactate m...	16272401	1
Caelo	Ethacridine lactate m...	17247606	1
Caelo	Ethacridine lactate m...	16272406	1
Caelo	Ethacridine lactate m...	162734	1
Caelo	Ethacridine lactate m...	15/214/R	1
Caelo	Ethacridine lactate m...	16273411	9
Caelo	Ethacridine lactate m...	17332314	3
Caelo	Ethacridine lactate m...	17332304	8
Caelo	Ethacridine lactate m...	7247805	1
Caelo	Ethacridine lactate m...	17332312	1
Caelo	Ethacridine lactate m...	141203404	1
Caelo	Ethacridine lactate m...	17247605	6
Caelo	Ethacridine lactate m...	1627412	1
Caelo	Ethacridine lactate m...	16272410	1
Caelo	Ethacridine lactate m...	16272408	3
Caelo	Ethacridine lactate m...	17332303	2
Caelo	Ethacridine lactate m...	12268407	3
Caelo	Eucalyptus leaf, cut	12225206	1
Caelo	Eucalyptus leaf, cut	122250204	1
Caelo	Eucalyptus leaf, cut	12225202	1
Caelo	Eucalyptus leaf, cut	12225203	3
Caelo	Eucalyptus leaf, cut	11307708	2
Caelo	Eucalyptus leaf, cut	9346811	1
Caelo	Euphrasy herb, cut	13180001	5
Caelo	Euphrasy herb, cut	8011419	1
Caelo	Euphrasy herb, cut	312052	1
Caelo	Euphrasy herb, cut	1213202	1
Caelo	Euphrasy herb, cut	7051302	1
Caelo	Euphrasy herb, cut	12103203	10
Caelo	Euphrasy herb, cut	13180004	5
Caelo	Euphrasy herb, cut	11030812	4
Caelo	Euphrasy herb, cut	12103202	9
Caelo	Euphrasy herb, cut	12103201	1
Caelo	Euphrasy herb, cut	29011302	1
Caelo	Euphrasy herb, cut	11030810	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Euphrasy herb, cut	13180005	2
Caelo	Euphrasy herb, cut	112211643	1
Caelo	Euphrasy herb, cut	11030811	3
Caelo	Euphrasy herb, cut	13180002	5
Caelo	Euphrasy herb, cut	30101307	1
Caelo	Euphrasy herb, cut	14019102	2
Caelo	Euphrasy herb, cut	14019107	1
Caelo	Euphrasy herb, cut	13180008	1
Caelo	European golden rod, ...	12343102	1
Caelo	European golden rod, ...	12343103-3466312121	1
Caelo	European golden rod, ...	13386908	1
Caelo	Fampridine	16276002	3
Caelo	Fampridine	16276001	1
Caelo	Fampridine	17217003	2
Caelo	Fampridine	24111801	1
Caelo	Fampridine	15350002	1
Caelo	Fampridine	15350001	1
Caelo	Fenugreek	12252602	2
Caelo	Fenugreek	1131206	1
Caelo	Fenugreek	13013805	3
Caelo	Fenugreek	11313203	1
Caelo	Fenugreek	12252603	3
Caelo	Fenugreek	12252601	2
Caelo	Fenugreek	13013803	6
Caelo	Fenugreek	13013801	1
Caelo	Fenugreek, powdered	12149605	1
Caelo	Fenugreek, powdered	7264a110201	2
Caelo	Fenugreek, powdered	12149607	1
Caelo	Fenugreek, powdered	13398302	1
Caelo	Fenugreek, powdered	14377201	2
Caelo	Fenugreek, powdered	14057801	1
Caelo	Fenugreek, powdered	14057807	1
Caelo	Fenugreek, powdered	15308703	2
Caelo	Feverfew herb, cut	12380204	1
Caelo	Feverfew herb, cut	12380203	5
Caelo	Finasteride (form III)	1809019-01	1
Caelo	Flax weed herb, cut	13007507	7
Caelo	Flax weed herb, cut	14005303	1
Caelo	Flax weed herb, cut	14005305	3
Caelo	Flax weed herb, cut	11104505	1
Caelo	Flax weed herb, cut	13007502	2
Caelo	Flax weed herb, cut	14197204	1
Caelo	Flax weed herb, cut	11104507	8
Caelo	Flax weed herb, cut	13007505	1
Caelo	Flax weed herb, cut	16292801	1
Caelo	Fluocinolonacetomid, ...	2804013V-03	1
Caelo	Fluocinolonacetomid, ...	16228108	2
Caelo	Fluocinolonacetomid, ...	12035003	1
Caelo	Frangula bark, cut	11254001	2
Caelo	Frangula bark, cut	11340801	1
Caelo	Frangula bark, cut	12185901	10
Caelo	Frangula bark, cut	12186002	15
Caelo	Frangula bark, cut	14163102	1
Caelo	Frangula bark, cut	12186001	1
Caelo	Frangula bark, cut	10061301	1
Caelo	Frangula bark, cut	6011402	1
Caelo	Frangula bark, cut	13203008	1
Caelo	Frangula bark, cut	12186005	1
Caelo	Frangula bark, cut	121860006	1
Caelo	Frangula bark, cut	12186006	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Frangula bark, cut	1286003	1
Caelo	Frangula bark, cut	12186004	11
Caelo	Frangula bark, cut	eztetdt	1
Caelo	Frangula bark, cut	12186003	1
Caelo	Frankincense	15140810	1
Caelo	Frankincense	15140805	2
Caelo	Frankincense	10349210	1
Caelo	Frankincense	12120005	5
Caelo	Frankincense	12120006	3
Caelo	Frankincense	15140809	1
Caelo	Frankincense	12120008	2
Caelo	Frankincense	12120012-3421310191	1
Caelo	Frankincense	10349213	1
Caelo	Frankincense	13112801	5
Caelo	Frankincense	13112803	1
Caelo	Frankincense	4614M-02030	1
Caelo	Frankincense	14107803	2
Caelo	Frankincense	13112809	2
Caelo	Frankincense	15140804	1
Caelo	Frankincense	13112807	1
Caelo	Frankincense	13112803-3494411030	1
Caelo	Fructose	15271214	2
Caelo	Fructose	18062005	1
Caelo	Fructose	18082005	13
Caelo	Fructose	18291605	1
Caelo	Fructose	13098913	5
Caelo	Fructose	13098922	1
Caelo	Fructose	10051401	1
Caelo	Fructose	13098915	1
Caelo	Fructose	13098917	3
Caelo	Fructose	13098909	1
Caelo	Fructose	31031401	1
Caelo	Fructose	13098911	1
Caelo	Fructose	13098914	2
Caelo	Fructose	1114A03394	1
Caelo	Fructose	31101810	1
Caelo	Fructose	12132320	1
Caelo	Fructose	13098918	1
Caelo	Fructose	15084312	2
Caelo	Fructose	17346411	2
Caelo	Fructose	18082002	4
Caelo	Fructose	1509037-01	12
Caelo	Fructose	14096201	2
Caelo	Fructose	15271205	5
Caelo	Fructose	14305206	1
Caelo	Fructose	15084311	1
Caelo	Fructose	15271213	6
Caelo	Fructose	15271208	1
Caelo	Fructose	15271201	1
Caelo	Fructose	15271209	3
Caelo	Fructose	16066103	3
Caelo	Fructose	15271215	3
Caelo	Fructose	16066106	9
Caelo	Fructose	16066107	7
Caelo	Fructose	23021708	1
Caelo	Fructose	17346405	1
Caelo	Fructose	17007701	5
Caelo	Fructose	16066101	2
Caelo	Fructose	170104901	1
Caelo	Fructose	17007708	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Fructose	5081701	1
Caelo	Fructose	17007707	2
Caelo	Fructose	1701034-01	1
Caelo	Fructose	2117E-03394	1
Caelo	Fructose	13111705	1
Caelo	Fructose	17143902	6
Caelo	Fructose	13041802	1
Caelo	Fructose	17143903	1
Caelo	Fructose	9061804	1
Caelo	Fructose	18061801	1
Caelo	Fructose	17346404	16
Caelo	Fructose	17007702	1
Caelo	Fructose	20081402	1
Caelo	Fructose	152711201	1
Caelo	Fructose	15084308	2
Caelo	Fructose	14096207	1
Caelo	Fructose	18021501	1
Caelo	Fructose	9031504	1
Caelo	Fructose	14096223	5
Caelo	Fructose	14305203	2
Caelo	Fructose	121e29-n08	1
Caelo	Fructose	14096211	1
Caelo	Fructose	14096210	1
Caelo	Fructose	12132314	5
Caelo	Fructose	14096224	2
Caelo	Fructose	5061503	1
Caelo	Fructose	14305201	1
Caelo	Fructose	3313U-03394	1
Caelo	Fructose	7051304	1
Caelo	Fructose	12132310	1
Caelo	Fructose	1081501	1
Caelo	Fructose	13098904	1
Caelo	Fructose	7111401	1
Caelo	Fructose	14096205	3
Caelo	Fructose	14096208	3
Caelo	Fructose	10381011	1
Caelo	Fructose	11258611	1
Caelo	Fructose	10071201	1
Caelo	Fructose	11071201	1
Caelo	Fructose	14096219	4
Caelo	Fructose	11258627	3
Caelo	Fructose	11258628	2
Caelo	Fructose	12132305	5
Caelo	Fructose	14096203	1
Caelo	Fructose	130998922	1
Caelo	Fructose	11258614	1
Caelo	Fructose	12132309	2
Caelo	Fructose	1309213394	1
Caelo	Fructose	15084305	3
Caelo	Fructose	14305204	1
Caelo	Fructose	21051303	1
Caelo	Fructose	12132322	2
Caelo	Fructose	13098902	3
Caelo	Fructose	12132325	3
Caelo	Fructose	12132313	5
Caelo	Fructose	12132317	2
Caelo	Fructose	15084301	3
Caelo	Fructose	1503045-01	9
Caelo	Fructose	11I07-N06	1
Caelo	Fructose	12132301	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Fumitory, cut	12103404	2
Caelo	Fumitory, cut	12103405	7
Caelo	Fumitory, cut	11101408	1
Caelo	Fumitory, cut	12103408	5
Caelo	Fumitory, cut	13202902	11
Caelo	Fumitory, cut	72878477	1
Caelo	Fumitory, cut	1320902	1
Caelo	Fumitory, cut	13202905	2
Caelo	Fumitory, cut	12103407	6
Caelo	Fumitory, cut	7111201	1
Caelo	Fumitory, cut	8312901	1
Caelo	Fumitory, cut	11101413	3
Caelo	Fumitory, cut	11101411	3
Caelo	Fumitory, cut	11101410	1
Caelo	Fumitory, cut	11101406	1
Caelo	Furosemide	16B15-B06	1
Caelo	Galangal root	15128008	1
Caelo	Galangal root	12330306	2
Caelo	Galangal root	16196904	1
Caelo	Galangal root	12330305	1
Caelo	Galangal root	11375904	1
Caelo	Galangal root	12125004	2
Caelo	Galangal root	1215004	1
Caelo	Galangal root	12330303	1
Caelo	Galangal root	16196901	1
Caelo	Garden violet, cut	11124008	1
Caelo	Garden violet, cut	11124007-3284310662	1
Caelo	Garden violet, cut	13038701	4
Caelo	Garden violet, cut	11124007	1
Caelo	Garden violet, cut	1303872	1
Caelo	Garden violet, cut	13038703	2
Caelo	Garden violet, cut	13232107	1
Caelo	Garden violet, cut	10011403	1
Caelo	Gelatine, white	12354401	1
Caelo	Gelatine, white	11334020	1
Caelo	Gelatine, white	11019308	1
Caelo	Gelatine, white	4110I-01578	1
Caelo	Gelatine, white	14130901	2
Caelo	Gelatine, white	14130912	1
Caelo	Gelatine, white	15084801	1
Caelo	Gelatine, white	15084806	1
Caelo	Gelatine, white	15320502	2
Caelo	Gelatine, white	8031701	1
Caelo	Gelatine, white	1544191916	1
Caelo	Gelatine, white	8051811	1
Caelo	Gelatine, white	15441922	1
Caelo	Gelatine, white	12354507	3
Caelo	Gelatine, white	11334025	4
Caelo	Gentamicin sulfate	150320G	1
Caelo	Gentamicin sulfate	143535	1
Caelo	Gentamicin sulfate	20161102-4	1
Caelo	Gentamicin sulfate	15226502	1
Caelo	Gentamicin sulfate	14353505	3
Caelo	Gentamicin sulfate	14K28-B06-303636	1
Caelo	Gentamicin sulfate	15258803	12
Caelo	Gentamicin sulfate	15258802	16
Caelo	Gentamicin sulfate	4563125	1
Caelo	Gentamicin sulfate	3714I-01582	1
Caelo	Gentamicin sulfate	16119302	1
Caelo	Gentamicin sulfate	14353504	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Gentamicin sulfate	FG1503073	1
Caelo	Gentamicin sulfate	15226508	2
Caelo	Gentamicin sulfate	15226513	1
Caelo	Gentamicin sulfate	15226518	5
Caelo	Gentamicin sulfate	1528802	1
Caelo	Gentamicin sulfate	1604054-02	9
Caelo	Gentamicin sulfate	15226502	1
Caelo	Gentamicin sulfate	ch.b.15258802	1
Caelo	Gentamicin sulfate	15258810	13
Caelo	Gentamicin sulfate	15258808	7
Caelo	Gentamicin sulfate	11113304	1
Caelo	Gentamicin sulfate	15258809	1
Caelo	Gentamicin sulfate	18300104	5
Caelo	Gentamicin sulfate	17332419	1
Caelo	Gentamicin sulfate	17358616	8
Caelo	Gentamicin sulfate	17358610	3
Caelo	Gentamicin sulfate	18G06-351916	1
Caelo	Gentamicin sulfate	17358604	5
Caelo	Gentamicin sulfate	FG1603072	1
Caelo	Gentamicin sulfate	17358603	13
Caelo	Gentamicin sulfate	12011801	1
Caelo	Gentamicin sulfate	18041804	1
Caelo	Gentamicin sulfate	1711025_01,180414G	1
Caelo	Gentamicin sulfate	17099606-20180404Ma1	1
Caelo	Gentamicin sulfate	17099608	11
Caelo	Gentamicin sulfate	161003602	1
Caelo	Gentamicin sulfate	1709601	1
Caelo	Gentamicin sulfate	17099604	12
Caelo	Gentamicin sulfate	17099606	9
Caelo	Gentamicin sulfate	1610036-01	23
Caelo	Gentamicin sulfate	1260417	1
Caelo	Gentamicin sulfate	15258815	3
Caelo	Gentamicin sulfate	15258816	16
Caelo	Gentamicin sulfate	15258819	5
Caelo	Gentamicin sulfate	15258818	22
Caelo	Gentamicin sulfate	15226501	18
Caelo	Gentamicin sulfate	15258819,20170818Ma1	1
Caelo	Gentamicin sulfate	20170829/1	1
Caelo	Gentamicin sulfate	15258805	1
Caelo	Gentamicin sulfate	115258818	1
Caelo	Gentamicin sulfate	1528818	1
Caelo	Gentamicin sulfate	17099603	17
Caelo	Gentamicin sulfate	17099601	2
Caelo	Gentamicin sulfate	15128818	1
Caelo	Gentamicin sulfate	143533504	1
Caelo	Gentamicin sulfate	13205007	2
Caelo	Gentamicin sulfate	1506016-01	12
Caelo	Gentamicin sulfate	13I19-B05-286565	1
Caelo	Gentamicin sulfate	13189608	3
Caelo	Gentamicin sulfate	13189611	11
Caelo	Gentamicin sulfate	13D24-N05	1
Caelo	Gentamicin sulfate	312012	1
Caelo	Gentamicin sulfate	73402189	1
Caelo	Gentamicin sulfate	100017814	1
Caelo	Gentamicin sulfate	13189607	8
Caelo	Gentamicin sulfate	13189606	3
Caelo	Gentamicin sulfate	13189602	5
Caelo	Gentamicin sulfate	13189603	3
Caelo	Gentamicin sulfate	6082013C	1
Caelo	Gentamicin sulfate	12146612	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Gentamicin sulfate	12146604	6
Caelo	Gentamicin sulfate	12146605	4
Caelo	Gentamicin sulfate	12163209	1
Caelo	Gentamicin sulfate	12146602	11
Caelo	Gentamicin sulfate	12163207	3
Caelo	Gentamicin sulfate	12G04-N16	6
Caelo	Gentamicin sulfate	1710I-01582	1
Caelo	Gentamicin sulfate	12163206	1
Caelo	Gentamicin sulfate	3110806	1
Caelo	Gentamicin sulfate	12163204	4
Caelo	Gentamicin sulfate	13091207	1
Caelo	Gentamicin sulfate	11113308	4
Caelo	Gentamicin sulfate	31032014E	1
Caelo	Gentamicin sulfate	13189606	3
Caelo	Gentamicin sulfate	13189612	9
Caelo	Gentamicin sulfate	13189612	1
Caelo	Gentamicin sulfate	14353503	1
Caelo	Gentamicin sulfate	12146602	1
Caelo	Gentamicin sulfate	14353504	11
Caelo	Gentamicin sulfate	17081502	1
Caelo	Gentamicin sulfate	14K28-B06-303646	8
Caelo	Gentamicin sulfate	1435303	1
Caelo	Gentamicin sulfate	14353502	3
Caelo	Gentamicin sulfate	1305047-03(2ma1)	1
Caelo	Gentamicin sulfate	13205005	1
Caelo	Gentamicin sulfate	130205005	1
Caelo	Gentamicin sulfate	14353503	13
Caelo	Gentamicin sulfate	17358619	1
Caelo	Gentamicin sulfate	15226502	13
Caelo	Gentamicin sulfate	13205006	1
Caelo	Gentamicin sulfate	13205006	5
Caelo	Gentamicin sulfate	13205004	1
Caelo	Gentamicin sulfate	23012015E	1
Caelo	Gentamicin sulfate	13205005	15
Caelo	Gentamicin sulfate	7011505	1
Caelo	Gentamicin sulfate	12146604	1
Caelo	Gentamicin sulfate	13189612	4
Caelo	Gentamicin sulfate	13189611	1
Caelo	Gentamicin sulfate	13205003	4
Caelo	Gentamicin sulfate	13205003	4
Caelo	Gentamicin sulfate	13189614	2
Caelo	Gentamicin sulfate	13189613	1
Caelo	Gentamicin sulfate	14E08-B03-294841	2
Caelo	Gentamicin sulfate	11113303	2
Caelo	Gentamicin sulfate	12146611	11
Caelo	Gentian root, cut	6311B130305	2
Caelo	Gentian root, cut	12244607	1
Caelo	Gentian root, cut	13156102	1
Caelo	Gentian root, cut	12244608	1
Caelo	Gentian root, cut	10263404	1
Caelo	Gentian root, cut	10263411	1
Caelo	Gentian root, cut	Int.Ch:313091	1
Caelo	Gentian root, cut	19071304	1
Caelo	Gentian root, cut	30011301	1
Caelo	Gentian root, powdered	12007106	1
Caelo	Gentian root, powdered	15011403	1
Caelo	Gentian root, powdered	14215601	1
Caelo	Ginger, peeled	13100604	1
Caelo	Ginger, peeled	11325712	1
Caelo	Ginger, peeled	11325716	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ginger, peeled	12180604	1
Caelo	Ginger, peeled	13100602	2
Caelo	Ginger, peeled	13310604	4
Caelo	Ginger, peeled	13310602	2
Caelo	Ginger, peeled	13100605	3
Caelo	Ginger, peeled	12164911-3368319002	1
Caelo	Ginger, peeled	11325703	2
Caelo	Ginger, peeled	19121201	1
Caelo	Ginger, peeled	11325715	2
Caelo	Ginger, peeled	11325719	3
Caelo	Ginger, peeled, powde...	13092013B	1
Caelo	Ginger, peeled, powde...	12164911	1
Caelo	Ginger, peeled, powde...	12164901	1
Caelo	Ginkgo leaf, cut	13340906	1
Caelo	Ginkgo leaf, cut	12373604	1
Caelo	Ginkgo leaf, cut	13340904	1
Caelo	Ginkgo leaf, cut	15077407	1
Caelo	Ginkgo leaf, cut	12151709	2
Caelo	Ginkgo leaf, cut	13340907	1
Caelo	Ginseng root	12080218	1
Caelo	Ginseng root	13341001	1
Caelo	Glucose	1415U-03393	1
Caelo	Glucose	16101205	1
Caelo	Glucose	15148808	4
Caelo	Glucose	15148803	1
Caelo	Glucose	15148810	4
Caelo	Glucose	712213393	1
Caelo	Glucose	11358823	1
Caelo	Glucose	12141715	1
Caelo	Glucose	16101204	1
Caelo	Glucose	15148807	3
Caelo	Glucose	14319308	1
Caelo	Glucose	1001313393	1
Caelo	Glucose	15329008	1
Caelo	Glucose	15329006	7
Caelo	Glucose	15329004	4
Caelo	Glucose	15329005	6
Caelo	Glucose	15148808	1
Caelo	Glucose	12294505	3
Caelo	Glucose	15148813	1
Caelo	Glucose	12294509	4
Caelo	Glucose	12294510	3
Caelo	Glucose	12141717	1
Caelo	Glucose	402313393	2
Caelo	Glucose	13178701	1
Caelo	Glucose	1509500	3
Caelo	Glucose	15329001	6
Caelo	Glucose	13329512	1
Caelo	Glucose	13329503	1
Caelo	Glucose	12141706	3
Caelo	Glucose	15095005	2
Caelo	Glucose	208223393	1
Caelo	Glucose	14142303	1
Caelo	Glucose	14319304	6
Caelo	Glucose	14142308	1
Caelo	Glucose	1914I03393	1
Caelo	Glucose	14051401	1
Caelo	Glucose	14033908	2
Caelo	Glucose	6031505	1
Caelo	Glucose	14033906	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glucose	14033907	5
Caelo	Glucose	14033901	2
Caelo	Glucose	5081401	1
Caelo	Glucose	3214M-03393	1
Caelo	Glucose	14142301	1
Caelo	Glucose	14142314	3
Caelo	Glucose	23091401	1
Caelo	Glucose	13041501	1
Caelo	Glucose	14033907	1
Caelo	Glucose	14319309	3
Caelo	Glucose	25061204	1
Caelo	Glucose	106213393	1
Caelo	Glucose	15095003	2
Caelo	Glucose	14142310	5
Caelo	Glucose	10091503	1
Caelo	Glucose	14142304	1
Caelo	Glucose	15148809	1
Caelo	Glucose	15148802	6
Caelo	Glucose	16050504	1
Caelo	Glucose	13329503	2
Caelo	Glucose	14319310	2
Caelo	Glucose	15095004	3
Caelo	Glucose	4071505	1
Caelo	Glucose	15095001	8
Caelo	Glucose	1439304	1
Caelo	Glucose	15148801	5
Caelo	Glucose	16050501	13
Caelo	Glucose	12141718	3
Caelo	Glucose	16050503	19
Caelo	Glucose	1812913	1
Caelo	Glucose	13178702	1
Caelo	Glucose	18012917	1
Caelo	Glucose	18012915	10
Caelo	Glucose	2021902	1
Caelo	Glucose	18012913	2
Caelo	Glucose	15329002	2
Caelo	Glucose	3513I-03393	1
Caelo	Glucose	17111803	1
Caelo	Glucose	18012912	5
Caelo	Glucose	12294519	1
Caelo	Glucose	12294517	1
Caelo	Glucose	13329512	6
Caelo	Glucose	1114A03393	1
Caelo	Glucose	13178704	3
Caelo	Glucose	14033902	5
Caelo	Glucose	12294514	2
Caelo	Glucose	13178707	2
Caelo	Glucose	18012921	2
Caelo	Glucose	19038002	1
Caelo	Glucose	13329504	2
Caelo	Glucose	13329501	3
Caelo	Glucose	214U-03393	1
Caelo	Glucose	5013I3393	1
Caelo	Glucose	1317815	1
Caelo	Glucose	13178718	1
Caelo	Glucose	4713L-03393	1
Caelo	Glucose	14033908	1
Caelo	Glucose	14142313	3
Caelo	Glucose	13178716	2
Caelo	Glucose	13178712	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glucose	3121402	1
Caelo	Glucose	12294509-3425311150	1
Caelo	Glucose	4313Q3393	1
Caelo	Glucose	AKD14	1
Caelo	Glucose	13178715	5
Caelo	Glucose	13178711	2
Caelo	Glucose	1614U3393	1
Caelo	Glucose	2513M03393	1
Caelo	Glucose	18012914	2
Caelo	Glucose	18012911	2
Caelo	Glucose	16328811	1
Caelo	Glucose	15021801	1
Caelo	Glucose	16328812	8
Caelo	Glucose	28111714	1
Caelo	Glucose	16328803	7
Caelo	Glucose	26041303	1
Caelo	Glucose	12141714	2
Caelo	Glucose	18012909	3
Caelo	Glucose	18111702	1
Caelo	Glucose	16328805	12
Caelo	Glucose	7091705	1
Caelo	Glucose	551616157406	1
Caelo	Glucose	16328807	4
Caelo	Glucose	6071717	1
Caelo	Glucose	16328801	12
Caelo	Glucose	2051701	1
Caelo	Glucose	16157411	7
Caelo	Glucose	17041701	1
Caelo	Glucose	8031703	1
Caelo	Glucose	16050508	1
Caelo	Glucose	16050507	2
Caelo	Glucose	16157406	18
Caelo	Glucose	160308G	1
Caelo	Glucose	16157403	27
Caelo	Glucose	1650502	1
Caelo	Glucose	16050502	12
Caelo	Glucose	1061703	1
Caelo	Glucose	12041801	1
Caelo	Glucose	14319302	1
Caelo	Glucose	12294515	2
Caelo	Glucose	18012907	9
Caelo	Glucose	14051802	1
Caelo	Glucose	15051301	1
Caelo	Glucose	8091810	1
Caelo	Glucose	15081831	1
Caelo	Glucose	2013Q-03393	1
Caelo	Glucose monohydrate	2315Q-02209NP	1
Caelo	Glucose monohydrate	15413405	4
Caelo	Glucose monohydrate	2101402	1
Caelo	Glucose monohydrate	16049103	3
Caelo	Glucose monohydrate	5092	1
Caelo	Glucose monohydrate	14102418	9
Caelo	Glucose monohydrate	15180001	10
Caelo	Glucose monohydrate	14102416	5
Caelo	Glucose monohydrate	1518002	1
Caelo	Glucose monohydrate	15095606	1
Caelo	Glucose monohydrate	16049101	20
Caelo	Glucose monohydrate	15413303	1
Caelo	Glucose monohydrate	15413412	4
Caelo	Glucose monohydrate	15095603	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glucose monohydrate	14334807	4
Caelo	Glucose monohydrate	15180015	7
Caelo	Glucose monohydrate	14102424	1
Caelo	Glucose monohydrate	15180003	3
Caelo	Glucose monohydrate	14102412	9
Caelo	Glucose monohydrate	14102411	4
Caelo	Glucose monohydrate	16049104	18
Caelo	Glucose monohydrate	123456	1
Caelo	Glucose monohydrate	13157013	1
Caelo	Glucose monohydrate	14102401	4
Caelo	Glucose monohydrate	1604104	1
Caelo	Glucose monohydrate	18071402	1
Caelo	Glucose monohydrate	150302G	1
Caelo	Glucose monohydrate	15180002	15
Caelo	Glucose monohydrate	16049108	4
Caelo	Glucose monohydrate	16184802	2
Caelo	Glucose monohydrate	16049105	9
Caelo	Glucose monohydrate	15180005	3
Caelo	Glucose monohydrate	14102415	1
Caelo	Glucose monohydrate	15095612	12
Caelo	Glucose monohydrate	14334808	2
Caelo	Glucose monohydrate	200415PJK01	1
Caelo	Glucose monohydrate	14334811	2
Caelo	Glucose monohydrate	14334806	3
Caelo	Glucose monohydrate	15413302	3
Caelo	Glucose monohydrate	15413402	7
Caelo	Glucose monohydrate	28051507	1
Caelo	Glucose monohydrate	14334810	6
Caelo	Glucose monohydrate	15413304/20160309/6	1
Caelo	Glucose monohydrate	15413410	4
Caelo	Glucose monohydrate	613011-3p	1
Caelo	Glucose monohydrate	13157022	1
Caelo	Glucose monohydrate	11021509	1
Caelo	Glucose monohydrate	14I04-B08-	1
Caelo	Glucose monohydrate	14102420	6
Caelo	Glucose monohydrate	14334802	15
Caelo	Glucose monohydrate	15413409	23
Caelo	Glucose monohydrate	160101	1
Caelo	Glucose monohydrate	15095602	1
Caelo	Glucose monohydrate	16049102	8
Caelo	Glucose monohydrate	15413413	1
Caelo	Glucose monohydrate	15095605	3
Caelo	Glucose monohydrate	15180010	8
Caelo	Glucose monohydrate	14102421	4
Caelo	Glucose monohydrate	15413401	5
Caelo	Glucose monohydrate	15413304	28
Caelo	Glucose monohydrate	15095611	3
Caelo	Glucose monohydrate	15180016	2
Caelo	Glucose monohydrate	15180012	16
Caelo	Glucose monohydrate	14102427	1
Caelo	Glucose monohydrate	14102430	5
Caelo	Glucose monohydrate	14334805	2
Caelo	Glucose monohydrate	17005813	6
Caelo	Glucose monohydrate	12284231	1
Caelo	Glucose monohydrate	18030003	6
Caelo	Glucose monohydrate	1707004-01	54
Caelo	Glucose monohydrate	11082014	1
Caelo	Glucose monohydrate	11082020	1
Caelo	Glucose monohydrate	1134506	1
Caelo	Glucose monohydrate	11314516	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glucose monohydrate	11L19-N13	4
Caelo	Glucose monohydrate	23051801	1
Caelo	Glucose monohydrate	1801009	1
Caelo	Glucose monohydrate	17305003	12
Caelo	Glucose monohydrate	17136804	9
Caelo	Glucose monohydrate	11314528	6
Caelo	Glucose monohydrate	12G30-N001	1
Caelo	Glucose monohydrate	11314526	1
Caelo	Glucose monohydrate	12284202	6
Caelo	Glucose monohydrate	17305004	9
Caelo	Glucose monohydrate	170058025517	1
Caelo	Glucose monohydrate	15180009	1
Caelo	Glucose monohydrate	12284204	2
Caelo	Glucose monohydrate	170058025621	1
Caelo	Glucose monohydrate	1803006	1
Caelo	Glucose monohydrate	18214710	1
Caelo	Glucose monohydrate	18214708	2
Caelo	Glucose monohydrate	UIUJHJH	1
Caelo	Glucose monohydrate	181197	1
Caelo	Glucose monohydrate	18214703	10
Caelo	Glucose monohydrate	18214704	7
Caelo	Glucose monohydrate	1810907	1
Caelo	Glucose monohydrate	18101907	12
Caelo	Glucose monohydrate	18214701	10
Caelo	Glucose monohydrate	18101909	1
Caelo	Glucose monohydrate	27081801	1
Caelo	Glucose monohydrate	24081801	1
Caelo	Glucose monohydrate	181401903	1
Caelo	Glucose monohydrate	18101901	18
Caelo	Glucose monohydrate	18101902	3
Caelo	Glucose monohydrate	18030005	4
Caelo	Glucose monohydrate	26061802	1
Caelo	Glucose monohydrate	18030006	1
Caelo	Glucose monohydrate	5051401	1
Caelo	Glucose monohydrate	170100401	1
Caelo	Glucose monohydrate	17136801	15
Caelo	Glucose monohydrate	13157005	12
Caelo	Glucose monohydrate	17005801	3
Caelo	Glucose monohydrate	12284230	4
Caelo	Glucose monohydrate	16184807	10
Caelo	Glucose monohydrate	16184808	6
Caelo	Glucose monohydrate	13157008	1
Caelo	Glucose monohydrate	E144C	2
Caelo	Glucose monohydrate	13157007	7
Caelo	Glucose monohydrate	13157004	1
Caelo	Glucose monohydrate	13157011	9
Caelo	Glucose monohydrate	1618401	2
Caelo	Glucose monohydrate	13157012	1
Caelo	Glucose monohydrate	16184801	16
Caelo	Glucose monohydrate	16049107	6
Caelo	Glucose monohydrate	2316E-02209	1
Caelo	Glucose monohydrate	4616M-02209	1
Caelo	Glucose monohydrate	13157019	10
Caelo	Glucose monohydrate	13157002	2
Caelo	Glucose monohydrate	16184816	3
Caelo	Glucose monohydrate	17005804	9
Caelo	Glucose monohydrate	13H06-N08	20
Caelo	Glucose monohydrate	170323183	1
Caelo	Glucose monohydrate	17005812	3
Caelo	Glucose monohydrate	17005805	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glucose monohydrate	12284206	5
Caelo	Glucose monohydrate	219317005803	1
Caelo	Glucose monohydrate	17136802	12
Caelo	Glucose monohydrate	170058025312	1
Caelo	Glucose monohydrate	170058025458	1
Caelo	Glucose monohydrate	17005811	2
Caelo	Glucose monohydrate	12284219	8
Caelo	Glucose monohydrate	17005802	12
Caelo	Glucose monohydrate	12284229	3
Caelo	Glucose monohydrate	2017Q-02209	1
Caelo	Glucose monohydrate	17005803	18
Caelo	Glucose monohydrate	11314519	1
Caelo	Glucose monohydrate	12284222	3
Caelo	Glucose monohydrate	14334809	1
Caelo	Glycerol monostearate...	13156601	2
Caelo	Glycerol monostearate...	12105704	1
Caelo	Glycerol monostearate...	12105729	1
Caelo	Glycerol monostearate...	14359407	1
Caelo	Glycerol monostearate...	15445007	5
Caelo	Glycerol monostearate...	18044701	2
Caelo	Glycerol monostearate...	15445014	1
Caelo	Glycerol monostearate...	18044702	1
Caelo	Glycerol monostearate...	15445011	1
Caelo	Glycerol monostearate...	15445010	1
Caelo	Glycerol monostearate...	15280912	1
Caelo	Glycerol monostearate...	15280904	2
Caelo	Glycerol monostearate...	13221308	1
Caelo	Glycerol monostearate...	131156602	1
Caelo	Glycerol monostearate...	14169710	1
Caelo	Glycerol monostearate...	13221303	2
Caelo	Glycerol monostearate...	16233804	6
Caelo	Glycerol monostearate...	13018201	3
Caelo	Glycerol monostearate...	1301821	1
Caelo	Glycerol monostearate...	18068303	2
Caelo	Glycerol monostearate...	12001001	1
Caelo	Glycerol monostearate...	13109304	12
Caelo	Glycerol monostearate...	17167804	3
Caelo	Glycerol monostearate...	15286307	2
Caelo	Glycerol monostearate...	Char. Bez. 17167803	1
Caelo	Glycerol monostearate...	17167803	4
Caelo	Glycerol monostearate...	15009203	16
Caelo	Glycerol monostearate...	15286302	3
Caelo	Glycerol monostearate...	15286308	3
Caelo	Glycerol monostearate...	17061703	3
Caelo	Glycerol monostearate...	14012002	5
Caelo	Glycerol monostearate...	1420607	1
Caelo	Glycerol monostearate...	14308310	4
Caelo	Glycerol monostearate...	14240607	1
Caelo	Glycerol monostearate...	14240611	1
Caelo	Glycerol monostearate...	16030304	2
Caelo	Glycerol monostearate...	18328702	2
Caelo	Glycerol monostearate...	16164606	1
Caelo	Glycerol monostearate...	17258501	1
Caelo	Glycerol monostearate...	17258505	1
Caelo	Glycerol monostearate...	13109410	1
Caelo	Glycerol monostearate...	16164604	1
Caelo	Glycerol monostearate...	13109402	1
Caelo	Glycerol monostearate...	16030301	3
Caelo	Glycerol monostearate...	13109405	2
Caelo	Glycerol monostearate...	11079308	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Glycine	15262503	1
Caelo	Glycine	15262504	1
Caelo	Glycine	15262502	2
Caelo	Glycine	14147005	1
Caelo	Glycine	17026502	1
Caelo	Glycine	15262518	1
Caelo	Goat's rue, cut	12343305	3
Caelo	Goat's rue, cut	12343302	4
Caelo	Goat's rue, cut	12017907	3
Caelo	Goat's rue, cut	5122013e	1
Caelo	Goat's rue, cut	12017912	1
Caelo	Goat's rue, cut	12343303	6
Caelo	Goat's rue, cut	11175104	1
Caelo	Goat's rue, cut	12343308	1
Caelo	Goat's rue, cut	11307901	1
Caelo	Goat's rue, cut	12017913	2
Caelo	Goat's rue, cut	12017914	1
Caelo	Goat's rue, cut	11307903	1
Caelo	Griseofulvin	100910	1
Caelo	Griseofulvin	16F13-B06-326074	1
Caelo	Ground ivy herb, cut	11374615	1
Caelo	Ground ivy herb, cut	11374609	1
Caelo	Ground ivy herb, cut	11052314	1
Caelo	Ground ivy herb, cut	11374606	1
Caelo	Ground ivy herb, cut	11052318	2
Caelo	Ground ivy herb, cut	11374602	1
Caelo	Ground ivy herb, cut	11374611	2
Caelo	Ground ivy herb, cut	11052315	1
Caelo	Guaiacum wood, cut	11352008	1
Caelo	Guaiacum wood, cut	11352010	1
Caelo	Guaiacum wood, cut	11352001	1
Caelo	Guaiacum wood, cut	11352012	1
Caelo	Guaiacum wood, cut	4011M-01778	1
Caelo	Guaiacum wood, cut	8041401	1
Caelo	Guaiacum wood, cut	16136710	1
Caelo	Guaiacum wood, cut	14008203	1
Caelo	Guaiacum wood, cut	16136706	1
Caelo	Guaifenesin	12277307	1
Caelo	Guaifenesin	1212002-01	1
Caelo	Guaifenesin	12277311	1
Caelo	Guaifenesin	15F30-b06-317816	1
Caelo	Guaifenesin	16135507	1
Caelo	Guaifenesin	15302502	1
Caelo	Hamamelis bark, cut	13123701	1
Caelo	Hamamelis bark, cut	12361301	1
Caelo	Hamamelis bark, cut	13123703	1
Caelo	Hamamelis leaf, cut	12164804	1
Caelo	Hamamelis leaf, cut	13269905	1
Caelo	Hamamelis leaf, cut	14356103	1
Caelo	Hamamelis leaf, cut	12164802	1
Caelo	Hamamelis leaf, cut	13269904	1
Caelo	Hamamelis leaf, cut	12165805	1
Caelo	Hamamelis leaf, cut	10102915	1
Caelo	Hamamelis leaf, cut	13269902	2
Caelo	Hard fat in pastilles	27101703	1
Caelo	Hard fat in pastilles	17307602	6
Caelo	Hard fat in pastilles	14338704	4
Caelo	Hard fat in pastilles	17107306	2
Caelo	Hard fat in pastilles	16274203	10
Caelo	Hard fat in pastilles	14338705	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hard fat in pastilles	17107301	6
Caelo	Hard fat in pastilles	170921AS/16274208	1
Caelo	Hard fat in pastilles	16274208	5
Caelo	Hard fat in pastilles	14338707	3
Caelo	Hard fat in pastilles	8051810	1
Caelo	Hard fat in pastilles	14338703	4
Caelo	Hard fat in pastilles	17307601	13
Caelo	Hard fat in pastilles	16274205	3
Caelo	Hard fat in pastilles	16027806	1
Caelo	Hard fat in pastilles	16274201	9
Caelo	Hard fat in pastilles	16027802	7
Caelo	Hard fat in pastilles	14126218	2
Caelo	Hard fat in pastilles	16027801	5
Caelo	Hard fat in pastilles	15185007	2
Caelo	Hard fat in pastilles	1518505	1
Caelo	Hard fat in pastilles	15185009	8
Caelo	Hard fat in pastilles	15185005	13
Caelo	Hard fat in pastilles	14338706	1
Caelo	Hard fat in pastilles	15185003	11
Caelo	Hard fat in pastilles	15185002	16
Caelo	Hard fat in pastilles	14338708	1
Caelo	Hard fat in pastilles	8041710	1
Caelo	Hard fat in pastilles	17307604	1
Caelo	Hard fat in pastilles	14338711	2
Caelo	Hard fat in pastilles	12K01-N12	3
Caelo	Hard fat in pastilles	7101304	1
Caelo	Hard fat in pastilles	18147908	1
Caelo	Hard fat in pastilles	13243411	3
Caelo	Hard fat in pastilles	16021503	1
Caelo	Hard fat in pastilles	13243409	5
Caelo	Hard fat in pastilles	13243419	3
Caelo	Hard fat in pastilles	11109623	2
Caelo	Hard fat in pastilles	11285001	1
Caelo	Hard fat in pastilles	14126202	4
Caelo	Hard fat in pastilles	14126203	9
Caelo	Hard fat in pastilles	13243412	2
Caelo	Hard fat in pastilles	13243420	1
Caelo	Hard fat in pastilles	13243414	1
Caelo	Hard fat in pastilles	12158115	8
Caelo	Hard fat in pastilles	11109617	2
Caelo	Hard fat in pastilles	13243416	1
Caelo	Hard fat in pastilles	90915429	1
Caelo	Hard fat in pastilles	12158114	3
Caelo	Hard fat in pastilles	B:11109606	1
Caelo	Hard fat in pastilles	11345905	3
Caelo	Hard fat in pastilles	12111204	1
Caelo	Hard fat in pastilles	13243418	3
Caelo	Hard fat in pastilles	12158122	1
Caelo	Hard fat in pastilles	11345918	9
Caelo	Hard fat in pastilles	13243301	6
Caelo	Hard fat in pastilles	1412602	1
Caelo	Hard fat in pastilles	15111402	1
Caelo	Hard fat in pastilles	18147902	3
Caelo	Hard fat in pastilles	12158109	3
Caelo	Hard fat in pastilles	14126217	3
Caelo	Hard fat in pastilles	10158813	1
Caelo	Hard fat in pastilles	29061803	1
Caelo	Hard fat in pastilles	14126216	2
Caelo	Hard fat in pastilles	11345917	3
Caelo	Hard fat in pastilles	1135917	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hard fat in pastilles	18147901	2
Caelo	Hard fat in pastilles	12158119	1
Caelo	Hard fat in pastilles	18040401	2
Caelo	Hard fat in pastilles	12158104	5
Caelo	Hard fat in pastilles	24011302	1
Caelo	Hard fat in pastilles	17041401	1
Caelo	Hard fat in pastilles	14126204	2
Caelo	Hard fat in pastilles	11345904	1
Caelo	Hard fat in pastilles	14126205	1
Caelo	Hard fat in pastilles	20071501	1
Caelo	Hard fat in pastilles	14126210	3
Caelo	Hard fat W35	18147901	1
Caelo	Hard fat W35	16274201	1
Caelo	Hard fat W35	16027805	1
Caelo	Hard fat W35	16K08-B10	1
Caelo	Hawthorn flower	13010402	2
Caelo	Hawthorn flower	7121201	1
Caelo	Hawthorn flower	1351B130819	3
Caelo	Hawthorn flower	11068208	1
Caelo	Hawthorn flower	11068202	1
Caelo	Hawthorn flower	11068206	1
Caelo	Hawthorn flower	13253304	1
Caelo	Hawthorn flower	25101301	1
Caelo	Hawthorn flower	13010403	4
Caelo	Hawthorn flower	11068209	1
Caelo	Hawthorn flower	21051301	1
Caelo	Hawthorn flower	13010402	3
Caelo	Hawthorn flower	13253303	1
Caelo	Hawthorn fruit	11278202	1
Caelo	Hawthorn fruit	12326306	1
Caelo	Hawthorn fruit	13367203	1
Caelo	Hawthorn fruit	12326304	1
Caelo	Hawthorn fruit	13101402	2
Caelo	Hawthorn fruit	12326307	2
Caelo	Hawthorn fruit	12326308	1
Caelo	Hawthorn leaf with fl. ...	11068206	1
Caelo	Hawthorn leaf with fl. ...	IntCh: 313082	1
Caelo	Hawthorn leaf with fl. ...	52467476	1
Caelo	Hawthorn leaf with fl. ...	1225107	1
Caelo	Hawthorn leaf with fl. ...	12225107	5
Caelo	Hawthorn leaf with fl. ...	1221707	1
Caelo	Hawthorn leaf with fl. ...	16025802	1
Caelo	Hawthorn leaf with fl. ...	1518802	1
Caelo	Hawthorn leaf with fl. ...	15188802	1
Caelo	Hawthorn leaf with fl. ...	12225108	7
Caelo	Hawthorn leaf with fl. ...	13276703	3
Caelo	Hawthorn leaf with fl. ...	1222514	1
Caelo	Hawthorn leaf with fl. ...	12225207	1
Caelo	Hawthorn leaf with fl. ...	109076	1
Caelo	Hawthorn leaf with fl. ...	10273512	2
Caelo	Hawthorn leaf with fl. ...	11366001	1
Caelo	Hawthorn leaf with fl. ...	2225104	1
Caelo	Hawthorn leaf with fl. ...	10273515	6
Caelo	Hawthorn leaf with fl. ...	12225104	12
Caelo	Hay flower	3713F-01459	1
Caelo	Hay flower	13097302	9
Caelo	Hay flower	25101202	1
Caelo	Hay flower	14121301	1
Caelo	Hay flower	11099902	1
Caelo	Hay flower	13099205	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hay flower	3713E-01459	1
Caelo	Hay flower	13097301	1
Caelo	Hay flower	12317110	2
Caelo	Hay flower	11099912	1
Caelo	Hay flower	12317101	8
Caelo	Hay flower	12317107	13
Caelo	Hay flower	12371707	1
Caelo	Hay flower	12317108	3
Caelo	Hay flower	1231702	1
Caelo	Hay flower	12317102	3
Caelo	Hay flower	11099915	1
Caelo	Hay flower	11099914	1
Caelo	Hay flower	11099909	1
Caelo	Hay flower	11099916	1
Caelo	Hay flower	12317104	10
Caelo	Heather flower, rubbed	10270516	1
Caelo	Heather flower, rubbed	13014206	1
Caelo	Heather flower, rubbed	13014204	4
Caelo	Heather, cut	28061201	1
Caelo	Heather, cut	13367502	1
Caelo	Heather, cut	19111203	1
Caelo	Heather, cut	11011209	3
Caelo	Heather, cut	11102901	1
Caelo	Heather, cut	9291512	1
Caelo	Heavy alkaline bismut...	10013213	1
Caelo	Heavy alkaline bismut...	11052805	1
Caelo	Heavy alkaline bismut...	11052809	1
Caelo	Heavy alkaline bismut...	10F09-N13	1
Caelo	Heavy alkaline bismut...	12146903	4
Caelo	Heavy alkaline bismut...	12146902	1
Caelo	Heavy alkaline bismut...	10013208	1
Caelo	Heavy alkaline bismut...	4513M-01186	1
Caelo	Heavy alkaline bismut...	12146913	2
Caelo	Heavy alkaline bismut...	M468105	1
Caelo	Heavy alkaline bismut...	12146911	3
Caelo	Henna powder	13262803	2
Caelo	Henna powder	12122310	3
Caelo	Henna powder	12267311	1
Caelo	Henna powder	12122301	2
Caelo	Histamine dihydrochlo...	130500601	1
Caelo	Hop strobile	12008309	3
Caelo	Hop strobile	14068007	1
Caelo	Hop strobile	10240308	1
Caelo	Hop strobile	12221804	1
Caelo	Hop strobile	6091303	1
Caelo	Hop strobile	13123803	7
Caelo	Hop strobile	12221807	1
Caelo	Hop strobile	13123902	1
Caelo	Hop strobile	13187602	1
Caelo	Hop strobile	13123802	9
Caelo	Hop strobile	4111304	1
Caelo	Hop strobile	30121301	1
Caelo	Hop strobile	13123803	1
Caelo	Hop strobile	13123804	3
Caelo	Hop strobile	131223805	1
Caelo	Hop strobile	13123803-3505411011	1
Caelo	Hop strobile	4756769	1
Caelo	Hop strobile	13123804-3563412081	1
Caelo	Hop strobile	13123807	1
Caelo	Hop strobile	13123806	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hop strobile	13123802	3
Caelo	Hop strobile	1451B130107	6
Caelo	Hop strobile	12008310	3
Caelo	Hop strobile	12008307	1
Caelo	Hop strobile	12008313	2
Caelo	Hop strobile	12092001	2
Caelo	Hop strobile	12006305	1
Caelo	Hop strobile	15061211	1
Caelo	Hop strobile	12008305	1
Caelo	Hop strobile	12221805	3
Caelo	Hop strobile	12008302	1
Caelo	Hop strobile	1451b120817	1
Caelo	Hop strobile	5021306	1
Caelo	Hop strobile	12008317	6
Caelo	Hop strobile	12008315	1
Caelo	Hop strobile	12380101	6
Caelo	Hop strobile	12008317	1
Caelo	Hop strobile	12380103	2
Caelo	Hop strobile	12221802	2
Caelo	Hop strobile	12008312	1
Caelo	Hop strobile	12008310	3
Caelo	Horse chestnut leaf, ...	11063503	2
Caelo	Huang Qi	g040h1202222	1
Caelo	Hyaluronic acid sodiu...	14364706	2
Caelo	Hyaluronic acid sodiu...	16246703	1
Caelo	Hyaluronic acid sodiu...	15448803	1
Caelo	Hyaluronic acid sodiu...	15204903	3
Caelo	Hyaluronic acid sodiu...	14188702	1
Caelo	Hyaluronic acid sodiu...	1263303	1
Caelo	Hyaluronic acid sodiu...	15204901	1
Caelo	Hydrochlorothiazide	15231009	1
Caelo	Hydrochlorothiazide	16269301	4
Caelo	Hydrochlorothiazide	16F17-B07-327489	2
Caelo	Hydrochlorothiazide	16269304	1
Caelo	Hydrochlorothiazide	16F17B07327489	1
Caelo	Hydrochlorothiazide	16269307	4
Caelo	Hydrochlorothiazide	16269309	2
Caelo	Hydrochlorothiazide	16269311	1
Caelo	Hydrochlorothiazide	18017303	1
Caelo	Hydrochlorothiazide	15231008	3
Caelo	Hydrochlorothiazide	18264901	2
Caelo	Hydrochlorothiazide	18017305	3
Caelo	Hydrochlorothiazide	15231011	1
Caelo	Hydrochlorothiazide	15231005	1
Caelo	Hydrochlorothiazide	18017304	3
Caelo	Hydrochlorothiazide	15231003	5
Caelo	Hydrochlorothiazide	11101n02	1
Caelo	Hydrochlorothiazide	13008507	1
Caelo	Hydrochlorothiazide	14009701	21
Caelo	Hydrochlorothiazide	15231002	3
Caelo	Hydrochlorothiazide	13008509	5
Caelo	Hydrochlorothiazide	16269313	2
Caelo	Hydrocortisone	18091507	9
Caelo	Hydrocortisone	810911	1
Caelo	Hydrocortisone	17351202,180523H	1
Caelo	Hydrocortisone	21101303	1
Caelo	Hydrocortisone	18091503	11
Caelo	Hydrocortisone	13069506	9
Caelo	Hydrocortisone	13069505	4
Caelo	Hydrocortisone	13155602	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrocortisone	1306906	1
Caelo	Hydrocortisone	12167010-3437311151	1
Caelo	Hydrocortisone	13069509	4
Caelo	Hydrocortisone	12104309	3
Caelo	Hydrocortisone	1216707	1
Caelo	Hydrocortisone	18091525	6
Caelo	Hydrocortisone	18091519	7
Caelo	Hydrocortisone	18091506	7
Caelo	Hydrocortisone	17256605	14
Caelo	Hydrocortisone	1756602	1
Caelo	Hydrocortisone	18091520	6
Caelo	Hydrocortisone	17351209-20180830Ma3	1
Caelo	Hydrocortisone	17256607	3
Caelo	Hydrocortisone	18091502	2
Caelo	Hydrocortisone	18091505	11
Caelo	Hydrocortisone	18091501-H20181115Ma	3
Caelo	Hydrocortisone	17351205	32
Caelo	Hydrocortisone	13155603	12
Caelo	Hydrocortisone	39/AKD17	1
Caelo	Hydrocortisone	18091508	2
Caelo	Hydrocortisone	1109845	1
Caelo	Hydrocortisone	7101301	1
Caelo	Hydrocortisone	2013U-01719	1
Caelo	Hydrocortisone	12104305	3
Caelo	Hydrocortisone	12k08-n01	1
Caelo	Hydrocortisone	17526602	1
Caelo	Hydrocortisone	12167007	10
Caelo	Hydrocortisone	12167001	9
Caelo	Hydrocortisone	12104308	8
Caelo	Hydrocortisone	11021302	1
Caelo	Hydrocortisone	12119303	9
Caelo	Hydrocortisone	12104307	17
Caelo	Hydrocortisone	16255608NP	1
Caelo	Hydrocortisone	3912Q-01719	1
Caelo	Hydrocortisone	113I-01719	1
Caelo	Hydrocortisone	3110456	1
Caelo	Hydrocortisone	13069501	8
Caelo	Hydrocortisone	12119304	3
Caelo	Hydrocortisone	12104304	9
Caelo	Hydrocortisone	11260309	3
Caelo	Hydrocortisone	K05B20160909	1
Caelo	Hydrocortisone	12104301	1
Caelo	Hydrocortisone	11260308	3
Caelo	Hydrocortisone	11273704	1
Caelo	Hydrocortisone	17351203	10
Caelo	Hydrocortisone	17351209	6
Caelo	Hydrocortisone	13069502	5
Caelo	Hydrocortisone	120403701	4
Caelo	Hydrocortisone	17351206	8
Caelo	Hydrocortisone	10282208	1
Caelo	Hydrocortisone	18091501	12
Caelo	Hydrocortisone	13063501	1
Caelo	Hydrocortisone	1726603	1
Caelo	Hydrocortisone	12104306	8
Caelo	Hydrocortisone	17119407	19
Caelo	Hydrocortisone	1626-B04-330340	1
Caelo	Hydrocortisone	1518012	1
Caelo	Hydrocortisone	16095201	21
Caelo	Hydrocortisone	14301304	1
Caelo	Hydrocortisone	16095202	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrocortisone	15021106	3
Caelo	Hydrocortisone	15021105	5
Caelo	Hydrocortisone	14b04-B02-293422	1
Caelo	Hydrocortisone	1501034-01	11
Caelo	Hydrocortisone	14301311	9
Caelo	Hydrocortisone	140301305	1
Caelo	Hydrocortisone	14301312	5
Caelo	Hydrocortisone	14301310	14
Caelo	Hydrocortisone	14115104	1
Caelo	Hydrocortisone	16095209	5
Caelo	Hydrocortisone	16095203	8
Caelo	Hydrocortisone	16095210	17
Caelo	Hydrocortisone	170208HC, 16004612	1
Caelo	Hydrocortisone	15/25-B01-318486	1
Caelo	Hydrocortisone	14301306	7
Caelo	Hydrocortisone	14115103	2
Caelo	Hydrocortisone	13342504	1
Caelo	Hydrocortisone	130617H	1
Caelo	Hydrocortisone	16095207	9
Caelo	Hydrocortisone	16004608	8
Caelo	Hydrocortisone	16004603	9
Caelo	Hydrocortisone	1518013	1
Caelo	Hydrocortisone	260416v4	1
Caelo	Hydrocortisone	15181012	15
Caelo	Hydrocortisone	15278702	6
Caelo	Hydrocortisone	15278703	8
Caelo	Hydrocortisone	15181013	12
Caelo	Hydrocortisone	15216101	6
Caelo	Hydrocortisone	15181007	5
Caelo	Hydrocortisone	15181014	6
Caelo	Hydrocortisone	15278705	10
Caelo	Hydrocortisone	15181006	18
Caelo	Hydrocortisone	15181003	7
Caelo	Hydrocortisone	16004604	24
Caelo	Hydrocortisone	15021104	6
Caelo	Hydrocortisone	16004607	15
Caelo	Hydrocortisone	15278706	6
Caelo	Hydrocortisone	16004605	10
Caelo	Hydrocortisone	15021107	13
Caelo	Hydrocortisone	3081604	1
Caelo	Hydrocortisone	14301307	1
Caelo	Hydrocortisone	14301313	4
Caelo	Hydrocortisone	16004611	10
Caelo	Hydrocortisone	17256604	20
Caelo	Hydrocortisone	14301305	5
Caelo	Hydrocortisone	16272305	8
Caelo	Hydrocortisone	16272308	4
Caelo	Hydrocortisone	16255680	1
Caelo	Hydrocortisone	16255608	12
Caelo	Hydrocortisone	17090003	4
Caelo	Hydrocortisone	17119402	13
Caelo	Hydrocortisone	17119404	1
Caelo	Hydrocortisone	170921HA/17090003	1
Caelo	Hydrocortisone	17090005	8
Caelo	Hydrocortisone	15181012/20160309/10	1
Caelo	Hydrocortisone	17119405	18
Caelo	Hydrocortisone	1672308	1
Caelo	Hydrocortisone	1719402-20171106Ma2	1
Caelo	Hydrocortisone	17119406	8
Caelo	Hydrocortisone	17119410	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrocortisone	17119403	17
Caelo	Hydrocortisone	4817Q-01719	1
Caelo	Hydrocortisone	17256602	28
Caelo	Hydrocortisone	167272306	1
Caelo	Hydrocortisone	13155604	6
Caelo	Hydrocortisone	17256804	1
Caelo	Hydrocortisone	17229410	1
Caelo	Hydrocortisone	16255604	2
Caelo	Hydrocortisone	14115105	3
Caelo	Hydrocortisone	16255603	6
Caelo	Hydrocortisone	16255607	4
Caelo	Hydrocortisone	14115108	2
Caelo	Hydrocortisone	14115107	9
Caelo	Hydrocortisone	13341511	3
Caelo	Hydrocortisone	16095212	5
Caelo	Hydrocortisone	13341508	5
Caelo	Hydrocortisone	13341510	12
Caelo	Hydrocortisone	13341509	9
Caelo	Hydrocortisone	13341501	6
Caelo	Hydrocortisone	13155614	4
Caelo	Hydrocortisone	13155611	2
Caelo	Hydrocortisone	13155605	5
Caelo	Hydrocortisone	31032014D	1
Caelo	Hydrocortisone	13341502	11
Caelo	Hydrocortisone	13155601	7
Caelo	Hydrocortisone	16004612/2020317	1
Caelo	Hydrocortisone	16095213	9
Caelo	Hydrocortisone	1140317	1
Caelo	Hydrocortisone	16272306	8
Caelo	Hydrocortisone	16272304	9
Caelo	Hydrocortisone	16255606	1
Caelo	Hydrocortisone	1672305	1
Caelo	Hydrocortisone	16095213-20170704Ma7	1
Caelo	Hydrocortisone	15181008	1
Caelo	Hydrocortisone	18091528	2
Caelo	Hydrocortisone	18091524	1
Caelo	Hydrocortisone acetate	15061006	1
Caelo	Hydrocortisone acetate	15356903	21
Caelo	Hydrocortisone acetate	16004402	14
Caelo	Hydrocortisone acetate	15067004,20160504Ma2	1
Caelo	Hydrocortisone acetate	15067013	10
Caelo	Hydrocortisone acetate	1430140	1
Caelo	Hydrocortisone acetate	1536902	1
Caelo	Hydrocortisone acetate	16004403	9
Caelo	Hydrocortisone acetate	16004407	17
Caelo	Hydrocortisone acetate	16004405	11
Caelo	Hydrocortisone acetate	20160727Ma2	1
Caelo	Hydrocortisone acetate	16004406	14
Caelo	Hydrocortisone acetate	16004408	16
Caelo	Hydrocortisone acetate	1600406	1
Caelo	Hydrocortisone acetate	160629HA	1
Caelo	Hydrocortisone acetate	160044402	1
Caelo	Hydrocortisone acetate	16095404	26
Caelo	Hydrocortisone acetate	16095406	21
Caelo	Hydrocortisone acetate	16004405,20161215Ma7	1
Caelo	Hydrocortisone acetate	16004407,20161215Ma9	1
Caelo	Hydrocortisone acetate	15H24_B01-316480	1
Caelo	Hydrocortisone acetate	16004404	1
Caelo	Hydrocortisone acetate	16095407	26
Caelo	Hydrocortisone acetate	16095405	17

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrocortisone acetate	16095408	3
Caelo	Hydrocortisone acetate	160044003	1
Caelo	Hydrocortisone acetate	16238704	5
Caelo	Hydrocortisone acetate	15356901	16
Caelo	Hydrocortisone acetate	15356904	8
Caelo	Hydrocortisone acetate	15356902	15
Caelo	Hydrocortisone acetate	15216204	2
Caelo	Hydrocortisone acetate	13291701	1
Caelo	Hydrocortisone acetate	14043904	10
Caelo	Hydrocortisone acetate	13291711	7
Caelo	Hydrocortisone acetate	14043901	10
Caelo	Hydrocortisone acetate	14043902	10
Caelo	Hydrocortisone acetate	14144806	11
Caelo	Hydrocortisone acetate	14144805	19
Caelo	Hydrocortisone acetate	14144803	7
Caelo	Hydrocortisone acetate	14144807	5
Caelo	Hydrocortisone acetate	14301405	26
Caelo	Hydrocortisone acetate	14301402	8
Caelo	Hydrocortisone acetate	14144802	1
Caelo	Hydrocortisone acetate	16238708	23
Caelo	Hydrocortisone acetate	14301406	19
Caelo	Hydrocortisone acetate	14301401	9
Caelo	Hydrocortisone acetate	14301407	12
Caelo	Hydrocortisone acetate	14301404	6
Caelo	Hydrocortisone acetate	14301408	5
Caelo	Hydrocortisone acetate	15A07-B03-306545	2
Caelo	Hydrocortisone acetate	15067004	20
Caelo	Hydrocortisone acetate	15067003	13
Caelo	Hydrocortisone acetate	15067005	18
Caelo	Hydrocortisone acetate	15067010	2
Caelo	Hydrocortisone acetate	301115	1
Caelo	Hydrocortisone acetate	15067011	6
Caelo	Hydrocortisone acetate	15067006	10
Caelo	Hydrocortisone acetate	1501035	1
Caelo	Hydrocortisone acetate	16238706	24
Caelo	Hydrocortisone acetate	166004407	1
Caelo	Hydrocortisone acetate	16238710	14
Caelo	Hydrocortisone acetate	17358404	1
Caelo	Hydrocortisone acetate	18278010	5
Caelo	Hydrocortisone acetate	25(1805010-01)	1
Caelo	Hydrocortisone acetate	18278003	5
Caelo	Hydrocortisone acetate	17219702	14
Caelo	Hydrocortisone acetate	17219709	7
Caelo	Hydrocortisone acetate	181115181	1
Caelo	Hydrocortisone acetate	180219HA/17219707	1
Caelo	Hydrocortisone acetate	18278004	9
Caelo	Hydrocortisone acetate	Caelo	2
Caelo	Hydrocortisone acetate	18142809	5
Caelo	Hydrocortisone acetate	18142807	3
Caelo	Hydrocortisone acetate	18278005	7
Caelo	Hydrocortisone acetate	18142804	13
Caelo	Hydrocortisone acetate	17219703	18
Caelo	Hydrocortisone acetate	18142805-H20181115Ma	1
Caelo	Hydrocortisone acetate	17219713	1
Caelo	Hydrocortisone acetate	17358403	5
Caelo	Hydrocortisone acetate	18142806	11
Caelo	Hydrocortisone acetate	18142812	7
Caelo	Hydrocortisone acetate	18142811	4
Caelo	Hydrocortisone acetate	18142810	1
Caelo	Hydrocortisone acetate	18142808	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrocortisone acetate	18142805	15
Caelo	Hydrocortisone acetate	172197055662	1
Caelo	Hydrocortisone acetate	18142802	5
Caelo	Hydrocortisone acetate	172197085556	1
Caelo	Hydrocortisone acetate	1403903	1
Caelo	Hydrocortisone acetate	18278014	1
Caelo	Hydrocortisone acetate	17219708	18
Caelo	Hydrocortisone acetate	16238714	11
Caelo	Hydrocortisone acetate	16238705	3
Caelo	Hydrocortisone acetate	16238712	28
Caelo	Hydrocortisone acetate	16238707	11
Caelo	Hydrocortisone acetate	1623708	1
Caelo	Hydrocortisone acetate	16238709	7
Caelo	Hydrocortisone acetate	1604052	1
Caelo	Hydrocortisone acetate	17086203	14
Caelo	Hydrocortisone acetate	16238713	7
Caelo	Hydrocortisone acetate	17086205	14
Caelo	Hydrocortisone acetate	170804HA/1604052-01	1
Caelo	Hydrocortisone acetate	17086205-20170818Ma3	1
Caelo	Hydrocortisone acetate	17219705	24
Caelo	Hydrocortisone acetate	16238713-20170818Ma2	1
Caelo	Hydrocortisone acetate	17053185	1
Caelo	Hydrocortisone acetate	17086206	9
Caelo	Hydrocortisone acetate	17086212	1
Caelo	Hydrocortisone acetate	17086207	3
Caelo	Hydrocortisone acetate	1708620402.2020	1
Caelo	Hydrocortisone acetate	163163055402	1
Caelo	Hydrocortisone acetate	171106HA/17086204	1
Caelo	Hydrocortisone acetate	17219704	18
Caelo	Hydrocortisone acetate	17219706	20
Caelo	Hydrocortisone acetate	17219706,20171127Ma2	1
Caelo	Hydrocortisone acetate	16095410	5
Caelo	Hydrocortisone acetate	170856206	1
Caelo	Hydrocortisone acetate	test	1
Caelo	Hydrocortisone acetate	14043907	13
Caelo	Hydrocortisone acetate	18142803	10
Caelo	Hydrocortisone acetate	14043903	14
Caelo	Hydrocortisone acetate	11263902	1
Caelo	Hydrocortisone acetate	1213M-01718	1
Caelo	Hydrocortisone acetate	12330106	4
Caelo	Hydrocortisone acetate	17219714	1
Caelo	Hydrocortisone acetate	12330103	9
Caelo	Hydrocortisone acetate	1205034-01	28
Caelo	Hydrocortisone acetate	12330102	8
Caelo	Hydrocortisone acetate	2612M01718	1
Caelo	Hydrocortisone acetate	7051303	1
Caelo	Hydrocortisone acetate	12330107	10
Caelo	Hydrocortisone acetate		3
Caelo	Hydrocortisone acetate	12330110	2
Caelo	Hydrocortisone acetate	12330104	9
Caelo	Hydrocortisone acetate	12330105	2
Caelo	Hydrocortisone acetate	12330112	13
Caelo	Hydrocortisone acetate	12330108	14
Caelo	Hydrocortisone acetate	12182408	2
Caelo	Hydrocortisone acetate	12330116	10
Caelo	Hydrocortisone acetate	3110491	1
Caelo	Hydrocortisone acetate	18011304	1
Caelo	Hydrocortisone acetate	1406211718	1
Caelo	Hydrocortisone acetate	15061207	1
Caelo	Hydrocortisone acetate	12008405	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydrocortisone acetate	11263905	1
Caelo	Hydrocortisone acetate	11098604	1
Caelo	Hydrocortisone acetate	11098615	1
Caelo	Hydrocortisone acetate	11263913	2
Caelo	Hydrocortisone acetate	12182405	21
Caelo	Hydrocortisone acetate	11263916	1
Caelo	Hydrocortisone acetate	11263917	3
Caelo	Hydrocortisone acetate	12008402	2
Caelo	Hydrocortisone acetate	12182403	5
Caelo	Hydrocortisone acetate	10301904	1
Caelo	Hydrocortisone acetate	12183403	1
Caelo	Hydrocortisone acetate	12182404	21
Caelo	Hydrocortisone acetate	12182402	11
Caelo	Hydrocortisone acetate	11263915	2
Caelo	Hydrocortisone acetate	12330111	3
Caelo	Hydrocortisone acetate	11098612	1
Caelo	Hydrocortisone acetate	13155306	4
Caelo	Hydrocortisone acetate	13155303	8
Caelo	Hydrocortisone acetate	1391704	1
Caelo	Hydrocortisone acetate	13155308	8
Caelo	Hydrocortisone acetate	1391706	1
Caelo	Hydrocortisone acetate	13291705	4
Caelo	Hydrocortisone acetate	3110516	1
Caelo	Hydrocortisone acetate	12330113	1
Caelo	Hydrocortisone acetate	131553080	1
Caelo	Hydrocortisone acetate	13155304	5
Caelo	Hydrocortisone acetate	13155307	2
Caelo	Hydrocortisone acetate	12330111-3507411061	1
Caelo	Hydrocortisone acetate	27091307	1
Caelo	Hydrocortisone acetate	13291704	3
Caelo	Hydrocortisone acetate	13291709	4
Caelo	Hydrocortisone acetate	13155302	9
Caelo	Hydrocortisone acetate	12330114	1
Caelo	Hydrocortisone acetate	13291707	14
Caelo	Hydrocortisone acetate	13291708	5
Caelo	Hydrocortisone acetate	13291706	9
Caelo	Hydrocortisone acetate	25101201	1
Caelo	Hydroquinone	17221508	3
Caelo	Hydroquinone	18017902	6
Caelo	Hydroquinone	16051502	1
Caelo	Hydroquinone	14280104	2
Caelo	Hydroquinone	2313I-01717Nachprüf	1
Caelo	Hydroquinone	20101806	1
Caelo	Hydroquinone	14280104	2
Caelo	Hydroquinone	15057001	7
Caelo	Hydroquinone	16023113	8
Caelo	Hydroquinone	4017A-01717	1
Caelo	Hydroquinone	16023112, 170808HC	1
Caelo	Hydroquinone	15023105	1
Caelo	Hydroquinone	160223105	1
Caelo	Hydroquinone	16023112	11
Caelo	Hydroquinone	52840166	1
Caelo	Hydroquinone	15057001	1
Caelo	Hydroquinone	16023105	1
Caelo	Hydroquinone	14026904	1
Caelo	Hydroquinone	8816-10-583	1
Caelo	Hydroquinone	15057002	1
Caelo	Hydroquinone	15057010	12
Caelo	Hydroquinone	1410401	1
Caelo	Hydroquinone	15057004	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydroquinone	15057004	1
Caelo	Hydroquinone	16023105	12
Caelo	Hydroquinone	12288505	14
Caelo	Hydroquinone	17221503	1
Caelo	Hydroquinone	14280105	2
Caelo	Hydroquinone	10145101	1
Caelo	Hydroquinone	1263912	1
Caelo	Hydroquinone	20200712-3	1
Caelo	Hydroquinone	14121202	1
Caelo	Hydroquinone	12156705	6
Caelo	Hydroquinone	11163912	5
Caelo	Hydroquinone	10145103	1
Caelo	Hydroquinone	11163905	4
Caelo	Hydroquinone	11163907	5
Caelo	Hydroquinone	12156703	2
Caelo	Hydroquinone	1228503	1
Caelo	Hydroquinone	12288502	2
Caelo	Hydroquinone	13105402	4
Caelo	Hydroquinone	16023118	7
Caelo	Hydroquinone	11163912/Int312052	1
Caelo	Hydroquinone	10145113	2
Caelo	Hydroquinone	13105408	9
Caelo	Hydroquinone	18031403	1
Caelo	Hydroquinone	14026902	5
Caelo	Hydroquinone	14180401	3
Caelo	Hydroquinone	21102014A	1
Caelo	Hydroquinone	14180410	4
Caelo	Hydroquinone	14180406	3
Caelo	Hydroquinone	14281501	11
Caelo	Hydroquinone	14281507	4
Caelo	Hydroquinone	1163902	1
Caelo	Hydroquinone	16023116	17
Caelo	Hydroquinone	14281504	1
Caelo	Hydroquinone	13105405	8
Caelo	Hydroquinone	18230602	1
Caelo	Hydroquinone	160231135430	1
Caelo	Hydroxyethylcellulose...	11121707	1
Caelo	Hydroxyethylcellulose...	15086505	1
Caelo	Hydroxyethylcellulose...	1337807	1
Caelo	Hydroxyethylcellulose...	140601701	1
Caelo	Hydroxyethylcellulose...	15086506	1
Caelo	Hydroxyethylcellulose...	15086504	2
Caelo	Hydroxyethylcellulose...	15086503	1
Caelo	Hydroxyethylcellulose...	1504018-1	1
Caelo	Hydroxyethylcellulose...	220914PML02	1
Caelo	Hydroxyethylcellulose...	1337808	1
Caelo	Hydroxyethylcellulose...	15086507	1
Caelo	Hydroxyethylcellulose...	16100606	17
Caelo	Hydroxyethylcellulose...	1406017-01	29
Caelo	Hydroxyethylcellulose...	13378515	3
Caelo	Hydroxyethylcellulose...	10019406	2
Caelo	Hydroxyethylcellulose...	16100605	3
Caelo	Hydroxyethylcellulose...	213M-03301	1
Caelo	Hydroxyethylcellulose...	19051406	1
Caelo	Hydroxyethylcellulose...	13378508	3
Caelo	Hydroxyethylcellulose...	17319606	1
Caelo	Hydroxyethylcellulose...	31031402	1
Caelo	Hydroxyethylcellulose...	1315306	1
Caelo	Hydroxyethylcellulose...	140317	1
Caelo	Hydroxyethylcellulose...	16100612	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydroxyethylcellulose...	13378502	4
Caelo	Hydroxyethylcellulose...	13153511	8
Caelo	Hydroxyethylcellulose...	10019402	2
Caelo	Hydroxyethylcellulose...	290114PML01	1
Caelo	Hydroxyethylcellulose...	12369718	2
Caelo	Hydroxyethylcellulose...	13153505	11
Caelo	Hydroxyethylcellulose...	13378510	1
Caelo	Hydroxyethylcellulose...	133785	1
Caelo	Hydroxyethylcellulose...	21011502	1
Caelo	Hydroxyethylcellulose...	13378506	3
Caelo	Hydroxyethylcellulose...	13378510	1
Caelo	Hydroxyethylcellulose...	15086405	3
Caelo	Hydroxyethylcellulose...	15086405	5
Caelo	Hydroxyethylcellulose...	19076708	1
Caelo	Hydroxyethylcellulose...	17319625	1
Caelo	Hydroxyethylcellulose...	19076703	1
Caelo	Hydroxyethylcellulose...	17319620	1
Caelo	Hydroxyethylcellulose...	17319615	1
Caelo	Hydroxyethylcellulose...	18C05-B05	1
Caelo	Hydroxyethylcellulose...	17319608	1
Caelo	Hydroxyethylcellulose...	13378524	1
Caelo	Hydroxyethylcellulose...	17319621	1
Caelo	Hydroxyethylcellulose...	171101201	8
Caelo	Hydroxyethylcellulose...	17319609	1
Caelo	Hydroxyethylcellulose...	12274720	1
Caelo	Hydroxyethylcellulose...	17J05-B03	1
Caelo	Hydroxyethylcellulose...	17319613	4
Caelo	Hydroxyethylcellulose...	17319603	2
Caelo	Hydroxyethylcellulose...	13378518	3
Caelo	Hydroxyethylcellulose...	14256002	2
Caelo	Hydroxyethylcellulose...	13378516	2
Caelo	Hydroxyethylcellulose...	15086403	3
Caelo	Hydroxyethylcellulose...	12011903	1
Caelo	Hydroxyethylcellulose...	13153511	1
Caelo	Hydroxyethylcellulose...	13041502	1
Caelo	Hydroxyethylcellulose...	15086407	6
Caelo	Hydroxyethylcellulose...	15086408	1
Caelo	Hydroxyethylcellulose...	1309020	1
Caelo	Hydroxyethylcellulose...	18113101	1
Caelo	Hydroxyethylcellulose...	17Jo5-B03-343927	1
Caelo	Hydroxyethylcellulose...	1731960608	1
Caelo	Hydroxyethylcellulose...	AKD10	1
Caelo	Hydroxyethylcellulose...	13378507	11
Caelo	Hydroxyethylcellulose...	12369717	1
Caelo	Hydroxyethylcellulose...	23061809	1
Caelo	Hydroxyethylcellulose...	16100615	4
Caelo	Hydroxyethylcellulose...	11266316	1
Caelo	Hydroxyethylcellulose...	16100620	3
Caelo	Hydroxyethylcellulose...	16100617	11
Caelo	Hydroxyethylcellulose...	1606030-02	28
Caelo	Hydroxyethylcellulose...	16H05325803	1
Caelo	Hydroxyethylcellulose...	16100610	3
Caelo	Hydroxyethylcellulose...	12156901	2
Caelo	Hydroxyethylcellulose...	16100622	4
Caelo	Hydroxyethylcellulose...	16100625	13
Caelo	Hydroxyethylcellulose...	19041708	1
Caelo	Hydroxyethylcellulose...	12156902	1
Caelo	Hydroxyethylcellulose...	21071802	1
Caelo	Hydroxyethylcellulose...	11266305	1
Caelo	Hydroxyethylcellulose...	13153506	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydroxyethylcellulose...	18022401	2
Caelo	Hydroxyethylcellulose...	700243-B04-331244	1
Caelo	Hydroxyethylcellulose...	12031803	1
Caelo	Hydroxyethylcellulose...	16278702	2
Caelo	Hydroxyethylcellulose...	8639291	1
Caelo	Hydroxyethylcellulose...	16278705	14
Caelo	Hydroxyethylcellulose...	19081705	1
Caelo	Hydroxyethylcellulose...	27101701	1
Caelo	Hydroxyethylcellulose...	16100390108639291	1
Caelo	Hydroxyethylcellulose...	16278714	7
Caelo	Hydroxyethylcellulose...	16278718	4
Caelo	Hydroxyethylcellulose...	16278715	2
Caelo	Hydroxyethylcellulose...	16100626	7
Caelo	Hydroxyethylcellulose...	12091819	1
Caelo	Hydroxyethylcellulose...	23041804	1
Caelo	Hydroxyethylcellulose...	12156907	9
Caelo	Hydroxyethylcellulose...	13153502	1
Caelo	Hydroxyethylcellulose...	26071301	1
Caelo	Hydroxyethylcellulose...	1110024-01	3
Caelo	Hydroxyethylcellulose...	13153501	10
Caelo	Hydroxyethylcellulose...	1304008-01	27
Caelo	Hydroxyethylcellulose...	12369713	2
Caelo	Hydroxyethylcellulose...	12249216	5
Caelo	Hydroxyethylcellulose...	122492014	1
Caelo	Hydroxyethylcellulose...	12369704	10
Caelo	Hydroxyethylcellulose...	8051305	1
Caelo	Hydroxyethylcellulose...	10019408	2
Caelo	Hydroxyethylcellulose...	20031305	1
Caelo	Hydroxyethylcellulose...	12249202	3
Caelo	Hydroxyethylcellulose...	25071301	1
Caelo	Hydroxyethylcellulose...	5012I-03301	1
Caelo	Hydroxyethylcellulose...	1113E-03301	1
Caelo	Hydroxyethylcellulose...	213I-03301	1
Caelo	Hydroxyethylcellulose...	1606030-01	43
Caelo	Hydroxyethylcellulose...	11186221	1
Caelo	Hydroxyethylcellulose...	21111812	1
Caelo	Hydroxyethylcellulose...	4512I-03301	1
Caelo	Hydroxyethylcellulose...	16100612	1
Caelo	Hydroxyethylcellulose...	16100607	1
Caelo	Hydroxyethylcellulose...	19011302	1
Caelo	Hydroxyethylcellulose...	16100613	2
Caelo	Hydroxyethylcellulose...	17319602	1
Caelo	Hydroxyethylcellulose...	12249205	5
Caelo	Hydroxyethylcellulose...	16100604	2
Caelo	Hydroxyethylcellulose...	16100606	1
Caelo	Hydroxyethylcellulose...	16252003	2
Caelo	Hydroxyethylcellulose...	15199510	1
Caelo	Hydroxyethylcellulose...	16252006	7
Caelo	Hydroxyethylcellulose...	16012901	1
Caelo	Hydroxyethylcellulose...	11304620	4
Caelo	Hydroxyethylcellulose...	14095901	1
Caelo	Hydroxyethylcellulose...	15199501	9
Caelo	Hydroxyethylcellulose...	14352510	2
Caelo	Hydroxyethylcellulose...	4091713	1
Caelo	Hydroxyethylcellulose...	13182702	6
Caelo	Hydroxyethylcellulose...	14352510	2
Caelo	Hydroxyethylcellulose...	13182714	1
Caelo	Hydroxyethylcellulose...	1309020-01	1
Caelo	Hydroxyethylcellulose...	14352508	1
Caelo	Hydroxyethylcellulose...	12274710	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hydroxyethylcellulose...	14095905	1
Caelo	Hydroxyethylcellulose...	18022502	2
Caelo	Hydroxyethylcellulose...	16252008	2
Caelo	Hydroxyethylcellulose...	23011901	1
Caelo	Hydroxyethylcellulose...	28111805	1
Caelo	Hydroxyethylcellulose...	28111804	1
Caelo	Hydroxyethylcellulose...	18022507	2
Caelo	Hydroxyethylcellulose...	16091711	1
Caelo	Hydroxyethylcellulose...	17101801	1
Caelo	Hydroxyethylcellulose...	14095911	4
Caelo	Hydroxyethylcellulose...	18022503	1
Caelo	Hydroxyethylcellulose...	9061802	1
Caelo	Hydroxyethylcellulose...	9041802	1
Caelo	Hydroxyethylcellulose...	17184503	8
Caelo	Hydroxyethylcellulose...	5121706	1
Caelo	Hydroxyethylcellulose...	18113101	3
Caelo	Hydroxypropylcellulos...	15345001	1
Caelo	Hydroxypropylcellulos...	15129707	1
Caelo	Hydroxypropylcellulos...	15345008	1
Caelo	Hydroxypropylcellulos...	13031309	1
Caelo	Hydroxypropylcellulos...	13031306	1
Caelo	Hydroxypropylcellulos...	14189801	1
Caelo	Hyoscine hydrobromide	16206208	1
Caelo	Hyoscine hydrobromide	15054511	1
Caelo	Hyoscine hydrobromide	16206204	1
Caelo	Hyoscine hydrobromide	15054505	1
Caelo	Hyoscine hydrobromide	14278801	26
Caelo	Hyoscine hydrobromide	18024807	1
Caelo	Hyoscine hydrobromide	18024806	3
Caelo	Hyoscine hydrobromide	17230303	1
Caelo	Hyoscine hydrobromide	15358901	1
Caelo	Hypromellose 4000	13d03-n03	1
Caelo	Hypromellose 4000	13155004	4
Caelo	Hypromellose 4000	12239502	1
Caelo	Hypromellose 4000	13100304	1
Caelo	Hypromellose 4000	12215505	1
Caelo	Hypromellose 4000	12215504	2
Caelo	Hypromellose 4000	12098702	1
Caelo	Hypromellose 4000	16224408	1
Caelo	Hypromellose 4000	13155002	1
Caelo	Hypromellose 4000	13100304	1
Caelo	Hypromellose 4000	13155004	6
Caelo	Hypromellose 4000	14338802	1
Caelo	Hypromellose 4000	15B02-B05-308852	2
Caelo	Hypromellose 4000	16224401	6
Caelo	Hypromellose 4000	16224403	14
Caelo	Hypromellose 4000	171150085501	1
Caelo	Hypromellose 4000	16224406	3
Caelo	Hypromellose 4000	16224407	6
Caelo	Hypromellose 4000	16224409	3
Caelo	Hypromellose 4000	180215H/17115004	1
Caelo	Hypromellose 4000	16224404	1
Caelo	Hypromellose 4000	14338801	3
Caelo	Hypromellose 4000	18132101	8
Caelo	Hypromellose 4000	17115005	16
Caelo	Hypromellose 4000	18132109	3
Caelo	Hypromellose 4000	18132102	2
Caelo	Hypromellose 4000	18132105	6
Caelo	Hypromellose 4000	19069601	2
Caelo	Hypromellose 4000	18132111	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Hypromellose 4000	17132101	1
Caelo	Hypromellose 4000	5121818	1
Caelo	Hyssop herb, cut	13143902	1
Caelo	Hyssop herb, cut	11374304	1
Caelo	Ibuprofen	17227101	1
Caelo	Ibuprofen	17227102	1
Caelo	Ibuprofen	17351401	2
Caelo	Ibuprofen	1410025.02	1
Caelo	Iceland moss, cut	12187002	1
Caelo	Iceland moss, cut	11053916	1
Caelo	Iceland moss, cut	14304402	1
Caelo	Iceland moss, cut	11053911	1
Caelo	Iceland moss, cut	14304401	1
Caelo	Iceland moss, cut	12197002	1
Caelo	Iceland moss, cut	12187102	20
Caelo	Iceland moss, cut	15297603	1
Caelo	Iceland moss, cut	14304403	2
Caelo	Iceland moss, cut	14304407	2
Caelo	Iceland moss, cut	14304409	6
Caelo	Iceland moss, cut	15297601	1
Caelo	Iceland moss, cut	12187103	1
Caelo	Iceland moss, cut	11101601	1
Caelo	Iceland moss, cut	15297605	2
Caelo	Iceland moss, cut	11101607	8
Caelo	Iceland moss, cut	9052804	2
Caelo	Iceland moss, cut	4912E-01774	1
Caelo	Iceland moss, cut	12187108	4
Caelo	Iceland moss, cut	13350807	4
Caelo	Iceland moss, cut	13350805	4
Caelo	Iceland moss, cut	12382806	1
Caelo	Iceland moss, cut	13350803	4
Caelo	Iceland moss, cut	12187107	1
Caelo	Iceland moss, cut	12187109	4
Caelo	Iceland moss, cut	12187106	2
Caelo	Iceland moss, cut	1287102	1
Caelo	Inositol	L1609025-01	1
Caelo	Iodoform	14077702	7
Caelo	Iodoform	14077603	2
Caelo	Iodoform	12283304	1
Caelo	Iodoform	12283303	1
Caelo	Iodoform	14077601	1
Caelo	Iodoform	12283206	3
Caelo	Iodoform	12283201	2
Caelo	Iodoform	14230001	5
Caelo	Iodoform	14229507	2
Caelo	Iodoform	14077607	1
Caelo	Iodoform	12150910	4
Caelo	Iodoform	14229504	3
Caelo	Iodoform	14230002	2
Caelo	Iodoform	14229508	3
Caelo	Iodoform	17086004	3
Caelo	Iodoform	14230006	4
Caelo	Iodoform	14077606	2
Caelo	Iodoform	12283202	5
Caelo	Iodoform	21051305	1
Caelo	Iodoform	12150905	2
Caelo	Iodoform	49(18091911)	1
Caelo	Iodoform	48(18091905)	1
Caelo	Iodoform	47(18091905)	1
Caelo	Iodoform	18091905	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Iodoform	18249502	1
Caelo	Iodoform	18091804	3
Caelo	Iodoform	18091801	1
Caelo	Iodoform	12070102	2
Caelo	Iodoform	12070502	1
Caelo	Iodoform	12070605	1
Caelo	Iodoform	12070501	1
Caelo	Iodoform	12070602	1
Caelo	Iodoform	12150902	5
Caelo	Iodoform	12150903	1
Caelo	Iodoform	12070606	2
Caelo	Iodoform	14230007	3
Caelo	Iodoform	12070608	3
Caelo	Iodoform	12150908	4
Caelo	Iodoform	14230003	1
Caelo	Iodoform	17086001	1
Caelo	Iodoform	14077602	1
Caelo	Iodoform	16012205	2
Caelo	Iodoform	16012204	3
Caelo	Iodoform	16012206	6
Caelo	Iodoform	14230008	4
Caelo	Iodoform	16012305	5
Caelo	Iodoform	16012207	4
Caelo	Iodoform	16012210	5
Caelo	Iodoform	15276113	1
Caelo	Iodoform	18091805	8
Caelo	Iodoform	14230304	2
Caelo	Iodoform	16012309	4
Caelo	Iodoform	12070105	1
Caelo	Iodoform	16012306	6
Caelo	Iodoform	16012308	3
Caelo	Iodoform	17086008	3
Caelo	Iodoform	17086005	7
Caelo	Iodoform	18091802	3
Caelo	Iodoform	17086002	2
Caelo	Iodoform	17086012	1
Caelo	Iodoform	17086006	9
Caelo	Iodoform	16012302	5
Caelo	Iodoform	10320407	1
Caelo	Iodoform	16012202	2
Caelo	Iodoform	16012201	5
Caelo	Iodoform	14230302	2
Caelo	Iodoform	14230305	7
Caelo	Iodoform	46(18091805)	1
Caelo	Iodoform	14230308	4
Caelo	Iodoform	14230306	4
Caelo	Iodoform	15276102	5
Caelo	Iodoform	15276001	3
Caelo	Iodoform	15276002	4
Caelo	Iodoform	18091806	5
Caelo	Iodoform	1010316	1
Caelo	Iodoform	15276101	5
Caelo	Iodoform	15276104	4
Caelo	Iodoform	1527611	1
Caelo	Iodoform	15276109	1
Caelo	Iodoform	15276103	4
Caelo	Iodoform	1527104	1
Caelo	Iodoform	15276107	1
Caelo	Iodoform	15276111	5
Caelo	Iodoform	1576103	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Iodoform	15276108	2
Caelo	Iodoform	15276105	6
Caelo	Iron (III) sulfate hy...	10354412	1
Caelo	Iron (III) sulfate hy...	10354415	1
Caelo	Iron (III) sulfate hy...	12175906	1
Caelo	Iron (III) sulfate hy...	12175908	1
Caelo	Iron (III) sulfate hy...	10354404	1
Caelo	Iron (III) sulfate hy...	1512210007	1
Caelo	Isosorbide dinitrate ...	14D09B03294281	1
Caelo	Isosorbide dinitrate ...	14D09-B03-307018	2
Caelo	Isosorbide dinitrate ...	14D-09-B03-300234	1
Caelo	Isosorbide dinitrate ...	13/26-N18	1
Caelo	Isosorbide dinitrate ...	13D15-N02	10
Caelo	Isosorbide dinitrate ...	17C29-B01-348860	1
Caelo	Isosorbide dinitrate ...	15-10-b01-328132	2
Caelo	Isosorbide dinitrate ...	17C29-B01-348859	12
Caelo	Isosorbide dinitrate ...	701319-0003	1
Caelo	Isosorbide dinitrate ...	15-10b01-328132	1
Caelo	Isosorbide dinitrate ...	17c29-b01-347885	4
Caelo	Isosorbide dinitrate ...	15/10-B01-325474	4
Caelo	Isosorbide dinitrate ...	1510b1320782	1
Caelo	Ispaghula seed	12276501	2
Caelo	Ispaghula seed	413E-03369	1
Caelo	Ispaghula seed	12276527	1
Caelo	Ispaghula seed	3051303	1
Caelo	Ispaghula seed	23011505	1
Caelo	Ispaghula seed	14246602	3
Caelo	Ispaghula seed	142660	1
Caelo	Ispaghula seed	15121901	1
Caelo	Ispaghula seed	15081503	1
Caelo	Ispaghula seed	10011303	1
Caelo	Ispaghula seed	18061804	1
Caelo	Ispaghula seed	15249319	2
Caelo	Ispaghula seed	28111703	1
Caelo	Ispaghula seed	15249315	1
Caelo	Ispaghula seed	29061510	1
Caelo	Ispaghula seed	4114001	1
Caelo	Ispaghula seed	4121202	1
Caelo	Ispaghula seed	12252709	1
Caelo	Ispaghula seed	30101302	1
Caelo	Ispaghula seed	16121302	1
Caelo	Ispaghula seed	26011401	1
Caelo	Ispaghula seed	5091319	1
Caelo	Ispaghula seed	9011302	1
Caelo	Ispaghula seed	1227506	1
Caelo	Ispaghula seed	12276511	1
Caelo	Ispaghula seed	19061308	1
Caelo	Ispaghula seed	11185811	2
Caelo	Ispaghula seed	14031408	1
Caelo	Ispaghula seed	4121201	1
Caelo	Ispaghula seed	11061405	1
Caelo	Ispaghula seed	18071301	1
Caelo	Ispaghula seed	8041702	1
Caelo	Ivy leaf, cut	27111301	1
Caelo	Ivy leaf, cut	12225605	5
Caelo	Ivy leaf, cut	11041303	1
Caelo	Ivy leaf, cut	12225602	1
Caelo	Ivy leaf, cut	2031301	1
Caelo	Ivy leaf, cut	11101205	1
Caelo	Ivy leaf, cut	8061205	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ivy leaf, cut	20111211	1
Caelo	Ivy leaf, cut	23121305	1
Caelo	Ivy leaf, cut	4011301	1
Caelo	Ivy leaf, cut	11229505	5
Caelo	Ivy leaf, cut	15101203	1
Caelo	Ivy leaf, cut	12225603	2
Caelo	Ivy leaf, cut	10011401	1
Caelo	Ivy leaf, cut	23121303	1
Caelo	Jasmine flower	12143203	3
Caelo	Java tea, cut	28101305	1
Caelo	Java tea, cut	2221B121112	16
Caelo	Java tea, cut	20111207	1
Caelo	Java tea, cut	2221A110825	2
Caelo	Java tea, cut	11212901	1
Caelo	Java tea, cut	12225702	1
Caelo	Java tea, cut	12225704	3
Caelo	Java tea, cut	21061304	1
Caelo	Java tea, cut	31121203	1
Caelo	Java tea, cut	13258804	3
Caelo	Java tea, cut	13258802	2
Caelo	Juniper	12028610	1
Caelo	Juniper	12330505	2
Caelo	Juniper	2890a101201	2
Caelo	Juniper, pounded	13031302	1
Caelo	Juniper, pounded	14149502	1
Caelo	Juniper, pounded	13367302	1
Caelo	Kaolin, heavy	10176518	1
Caelo	Kaolin, heavy	11350125	4
Caelo	Kaolin, heavy	11350124	2
Caelo	Kaolin, heavy	11350107	1
Caelo	Kaolin, heavy	13041803	1
Caelo	Kaolin, heavy	13041803	1
Caelo	Kaolin, heavy	63033516	1
Caelo	Kaolin, heavy	13041835	1
Caelo	Kaolin, heavy	13041835	1
Caelo	Kaolin, heavy	82043269	1
Caelo	Kaolin, heavy	11041610	1
Caelo	Kelp, cut	14071502	1
Caelo	Kelp, cut	11351109	1
Caelo	Kelp, cut	280172	4
Caelo	Kelp, cut	11352209	1
Caelo	Ketoconazole	13179612	2
Caelo	Ketoconazole	13399404	2
Caelo	Ketoconazole	13399404	6
Caelo	Ketoconazole	13399408	5
Caelo	Ketoconazole	19096601	1
Caelo	Ketoconazole	12300406	1
Caelo	Ketoconazole	13399408	1
Caelo	Ketoconazole	53317046	1
Caelo	Ketoconazole	14350420	2
Caelo	Ketoconazole	18007020	1
Caelo	Ketoconazole	53317046	1
Caelo	Ketoconazole	13179601	1
Caelo	Ketoconazole	13179606	3
Caelo	Ketoconazole	18007011	1
Caelo	Ketoconazole	13399412	1
Caelo	Ketoconazole	14350434	1
Caelo	Ketoconazole	14350419	1
Caelo	Ketoconazole	14350436	2
Caelo	Ketoconazole	14350437	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ketoconazole	14350427	2
Caelo	Ketoconazole	14350422	3
Caelo	Ketoconazole	14350421	1
Caelo	Ketoconazole	12300407	2
Caelo	Ketoconazole	14350416	3
Caelo	Ketoconazole	11070022	1
Caelo	Ketoconazole	15347701	7
Caelo	Ketoconazole	14350442	3
Caelo	Ketoconazole	14k04-b06-309671	1
Caelo	Ketoconazole	18007003	2
Caelo	Ketoconazole	14350408	2
Caelo	Ketoconazole	11070001	1
Caelo	Ketoconazole	14350429	3
Caelo	Kidney vetch, coarsel...	13104301	2
Caelo	Knotgrass, cut	4312M-01674	1
Caelo	Knotgrass, cut	7102427	1
Caelo	Knotgrass, cut	14079304	1
Caelo	Knotgrass, cut	12368903	1
Caelo	Knotgrass, cut	12368906	3
Caelo	Knotgrass, cut	12295704	4
Caelo	Knotgrass, cut	12295702	3
Caelo	Knotgrass, cut	8111202	1
Caelo	Knotgrass, cut	4301A100126	2
Caelo	Knotgrass, cut	12368902	5
Caelo	Knotgrass, cut	12368904	1
Caelo	Lactose monohydrate	16027708	8
Caelo	Lactose monohydrate	s-100513-1	1
Caelo	Lactose monohydrate	143335114	1
Caelo	Lactose monohydrate	14333521	7
Caelo	Lactose monohydrate	15081701	1
Caelo	Lactose monohydrate	16091710	1
Caelo	Lactose monohydrate	1703019-01	9
Caelo	Lactose monohydrate	6101702	1
Caelo	Lactose monohydrate	17078307	4
Caelo	Lactose monohydrate	17078302	6
Caelo	Lactose monohydrate	17078306	7
Caelo	Lactose monohydrate	27101702	1
Caelo	Lactose monohydrate	14333518	5
Caelo	Lactose monohydrate	8071501	1
Caelo	Lactose monohydrate	23061504	1
Caelo	Lactose monohydrate	14333520	8
Caelo	Lactose monohydrate	1433510	1
Caelo	Lactose monohydrate	16231215	4
Caelo	Lactose monohydrate	1433502	1
Caelo	Lactose monohydrate	5061504	1
Caelo	Lactose monohydrate	1501017-01	16
Caelo	Lactose monohydrate	14333509	11
Caelo	Lactose monohydrate	23051501	1
Caelo	Lactose monohydrate	14333510	1
Caelo	Lactose monohydrate	17078303	22
Caelo	Lactose monohydrate	112192	2
Caelo	Lactose monohydrate	14333511	5
Caelo	Lactose monohydrate	17078312	9
Caelo	Lactose monohydrate	162312135465	1
Caelo	Lactose monohydrate	14333501	3
Caelo	Lactose monohydrate	14333513	5
Caelo	Lactose monohydrate	1703019-020	1
Caelo	Lactose monohydrate	14333514	1
Caelo	Lactose monohydrate	16231212	5
Caelo	Lactose monohydrate	15236217	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lactose monohydrate	15236220	4
Caelo	Lactose monohydrate	1433520	2
Caelo	Lactose monohydrate	15236213	2
Caelo	Lactose monohydrate	16231201	4
Caelo	Lactose monohydrate	15236209	7
Caelo	Lactose monohydrate	1526201	1
Caelo	Lactose monohydrate	15236203	2
Caelo	Lactose monohydrate	16231207	7
Caelo	Lactose monohydrate	1090216	1
Caelo	Lactose monohydrate	123456	1
Caelo	Lactose monohydrate	16231208	5
Caelo	Lactose monohydrate	15236214	8
Caelo	Lactose monohydrate	16231213	4
Caelo	Lactose monohydrate	5041702	1
Caelo	Lactose monohydrate	16231211	23
Caelo	Lactose monohydrate	1623120803815470	1
Caelo	Lactose monohydrate	16231206	6
Caelo	Lactose monohydrate	16321202	1
Caelo	Lactose monohydrate	15011701	1
Caelo	Lactose monohydrate	1433518	2
Caelo	Lactose monohydrate	14333528	2
Caelo	Lactose monohydrate	2051702	1
Caelo	Lactose monohydrate	16231214	5
Caelo	Lactose monohydrate	10071701	1
Caelo	Lactose monohydrate	15236201	21
Caelo	Lactose monohydrate	15236204	10
Caelo	Lactose monohydrate	1433521	1
Caelo	Lactose monohydrate	15236208	12
Caelo	Lactose monohydrate	15236215	4
Caelo	Lactose monohydrate	17321001	18
Caelo	Lactose monohydrate	17321004	9
Caelo	Lactose monohydrate	5011901	1
Caelo	Lactose monohydrate	17336013	5
Caelo	Lactose monohydrate	17336014	2
Caelo	Lactose monohydrate	18198001	5
Caelo	Lactose monohydrate	1733008	1
Caelo	Lactose monohydrate	3031501	1
Caelo	Lactose monohydrate	18198007	5
Caelo	Lactose monohydrate	18198002	7
Caelo	Lactose monohydrate	14025608	2
Caelo	Lactose monohydrate	14140309	2
Caelo	Lactose monohydrate	18198008	1
Caelo	Lactose monohydrate	18198014	1
Caelo	Lactose monohydrate	14025619	3
Caelo	Lactose monohydrate	2011901	1
Caelo	Lactose monohydrate	14140308	6
Caelo	Lactose monohydrate	7111402	1
Caelo	Lactose monohydrate	14140305	6
Caelo	Lactose monohydrate	23101404	1
Caelo	Lactose monohydrate	9051401	1
Caelo	Lactose monohydrate	14025620	13
Caelo	Lactose monohydrate	1402504	1
Caelo	Lactose monohydrate	14140302	4
Caelo	Lactose monohydrate	14140304	2
Caelo	Lactose monohydrate	3091405	1
Caelo	Lactose monohydrate	14025603	7
Caelo	Lactose monohydrate	14025616	6
Caelo	Lactose monohydrate	26071401	1
Caelo	Lactose monohydrate	140612	3
Caelo	Lactose monohydrate	3121406	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lactose monohydrate	14333508	1
Caelo	Lactose monohydrate	3788916	1
Caelo	Lactose monohydrate	14333503	1
Caelo	Lactose monohydrate	14025604	9
Caelo	Lactose monohydrate	13299102	15
Caelo	Lactose monohydrate	24011405	1
Caelo	Lactose monohydrate	8021401	1
Caelo	Lactose monohydrate	20021411	1
Caelo	Lactose monohydrate	13299107	5
Caelo	Lactose monohydrate	13299108	8
Caelo	Lactose monohydrate	7031402	1
Caelo	Lactose monohydrate	12399107	1
Caelo	Lactose monohydrate	13299103	1
Caelo	Lactose monohydrate	12041401	1
Caelo	Lactose monohydrate	200415HJK02	1
Caelo	Lactose monohydrate	11041801	1
Caelo	Lactose monohydrate	21111813	1
Caelo	Lactose monohydrate	17321002	3
Caelo	Lactose monohydrate	9061803	1
Caelo	Lactose monohydrate	17078303-20180619Ma2	1
Caelo	Lactose monohydrate	17336005	4
Caelo	Lactose monohydrate	26061801	1
Caelo	Lactose monohydrate	30031503	1
Caelo	Lactose monohydrate	17336002	7
Caelo	Lactose monohydrate	17336001	1
Caelo	Lactose monohydrate	21071801	1
Caelo	Lactose monohydrate	17336009	4
Caelo	Lactose monohydrate	17336007	20
Caelo	Lactose monohydrate	30081801	1
Caelo	Lactose monohydrate	17336008	6
Caelo	Lactose monohydrate	14333502	2
Caelo	Lactose monohydrate	17321005	1
Caelo	Lactose monohydrate	16231202	17
Caelo	Lactose monohydrate	14025601	4
Caelo	Lactose monohydrate	1071302	1
Caelo	Lactose monohydrate	11219216	2
Caelo	Lactose monohydrate	16027707	2
Caelo	Lactose monohydrate	150903002	2
Caelo	Lactose monohydrate	13104909-3422310182	1
Caelo	Lactose monohydrate	11079002	2
Caelo	Lactose monohydrate	13104909	7
Caelo	Lactose monohydrate	16027706	5
Caelo	Lactose monohydrate	11219208	3
Caelo	Lactose monohydrate	4313Q2210	1
Caelo	Lactose monohydrate	13104908	2
Caelo	Lactose monohydrate	123Test	1
Caelo	Lactose monohydrate	12021508	1
Caelo	Lactose monohydrate	160913PJK01	1
Caelo	Lactose monohydrate	13081302	1
Caelo	Lactose monohydrate	10195107	2
Caelo	Lactose monohydrate	12336814	1
Caelo	Lactose monohydrate	13104902	7
Caelo	Lactose monohydrate	152362	2
Caelo	Lactose monohydrate	12336820	2
Caelo	Lactose monohydrate	16027704	1
Caelo	Lactose monohydrate	123test	1
Caelo	Lactose monohydrate	16027702	6
Caelo	Lactose monohydrate	16027701	3
Caelo	Lactose monohydrate	1523614	1
Caelo	Lactose monohydrate	6091301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lactose monohydrate	12021518	5
Caelo	Lactose monohydrate	12021502	1
Caelo	Lactose monohydrate	30101201	1
Caelo	Lactose monohydrate	15236202	1
Caelo	Lactose monohydrate	12336816	2
Caelo	Lactose monohydrate	16027703	6
Caelo	Lactose monohydrate	13299101	9
Caelo	Lactose monohydrate	16231203	4
Caelo	Lactose monohydrate	12336815	14
Caelo	Lactose monohydrate	12139909	1
Caelo	Lactose monohydrate	13104912	1
Caelo	Lactose monohydrate	12139912	1
Caelo	Lactose monohydrate	12336801	8
Caelo	Lactose monohydrate		1
Caelo	Lactose monohydrate	1301C005	2
Caelo	Lactose monohydrate	12139920	6
Caelo	Lactose monohydrate	12139910	1
Caelo	Lactose monohydrate	13104909-3453311172	1
Caelo	Lactose monohydrate	12139901	2
Caelo	Lactose monohydrate	12139908	1
Caelo	Lactose monohydrate	12139914	7
Caelo	Lactose monohydrate	11014714	1
Caelo	Lactose monohydrate	17011301	1
Caelo	Lactose monohydrate	2901312210	1
Caelo	Lactose monohydrate	13104911	2
Caelo	Lactose monohydrate	1106022-02	1
Caelo	Lactose monohydrate	13104910	2
Caelo	Lactose monohydrate	15236221	6
Caelo	Ladys bedstraw herb,...	12369007	3
Caelo	Ladys bedstraw herb,...	15029906	1
Caelo	Ladys bedstraw herb,...	11265104	1
Caelo	Ladys bedstraw herb,...	12369001	3
Caelo	Ladys bedstraw herb,...	12369006	5
Caelo	Ladys bedstraw herb,...	2401311647	1
Caelo	Ladys bedstraw herb,...	7011301	1
Caelo	Ladys bedstraw herb,...	11265102	9
Caelo	Ladys bedstraw herb,...	11265106	3
Caelo	Ladys bedstraw herb,...	11265103	5
Caelo	Ladys bedstraw herb,...	24101201	1
Caelo	Ladys bedstraw herb,...	10375707	2
Caelo	Ladys bedstraw herb,...	14019202	1
Caelo	Ladys bedstraw herb,...	12369002	3
Caelo	Ladys bedstraw herb,...	22041303	1
Caelo	Ladys bedstraw herb,...	10375702	1
Caelo	Lapacho wood, cut	2x11154506	1
Caelo	Lapacho wood, cut	12239303	2
Caelo	Lapacho wood, cut	13298308	1
Caelo	Lapacho wood, cut	13298305	1
Caelo	Lapacho wood, cut	11154506	1
Caelo	Lapacho wood, cut	19061301	1
Caelo	Lapacho wood, cut	12239305	8
Caelo	Large waybread herb, ...	9311208	1
Caelo	Large waybread herb, ...	13232401	1
Caelo	Large waybread herb, ...	10349610	1
Caelo	Large waybread herb, ...	14125907	1
Caelo	Large waybread herb, ...	10349603	2
Caelo	Large waybread herb, ...	10349614	1
Caelo	Large waybread herb, ...	13232404	1
Caelo	Large waybread herb, ...	12374106	2
Caelo	Large waybread herb, ...	12374101	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Large waybread herb, ...	12374104	7
Caelo	Large waybread herb, ...	123741004	1
Caelo	Laurel leaf	12299101	1
Caelo	Laurel leaf	12299106	3
Caelo	Lavender flower	13185001-3415310152	1
Caelo	Lavender flower	13185001-3389311070	1
Caelo	Lavender flower	11327903	5
Caelo	Lavender flower	15061209	1
Caelo	Lavender flower	27091202	1
Caelo	Lavender flower	13185002	1
Caelo	Lavender flower	12031302	1
Caelo	Lavender flower	13185001	27
Caelo	Lavender flower	12088002	4
Caelo	Lavender flower	12088003	8
Caelo	Lavender flower	13185006	1
Caelo	Lavender flower	1318501	1
Caelo	Lavender flower	15061201	1
Caelo	Lavender flower	19071303	1
Caelo	Lavender flower	13185003	10
Caelo	Lavender flower	12088005	3
Caelo	Lavender flower	6091304	1
Caelo	Lavender flower	12088006	3
Caelo	Lavender flower	11327905	2
Caelo	Lavender flower	13336008	1
Caelo	Lavender flower	13159601	1
Caelo	Lavender flower	13387204	1
Caelo	Lavender flower	11054104	2
Caelo	Lavender flower	11054106	2
Caelo	Lavender flower	11054107	1
Caelo	Lemon grass, cut	11301501	1
Caelo	Lemon grass, cut	11301503	1
Caelo	Lemon grass, cut	13149701	4
Caelo	Lemon grass, cut	11301509	2
Caelo	Lemon grass, cut	13149703	1
Caelo	Lemon grass, cut	13341301	1
Caelo	Lemon peel, cut	12217401	1
Caelo	Lemon peel, cut	12217402	1
Caelo	Lemon peel, cut	12217404	1
Caelo	Lemon verbena leaf, cut	12137507	3
Caelo	Lemon verbena leaf, cut	10018208	1
Caelo	Lemon verbena leaf, cut	13037702	1
Caelo	Lemon verbena leaf, cut	12137505	1
Caelo	Lidocaine	12163810	1
Caelo	Lidocaine	18041812	1
Caelo	Lidocaine	16289503	6
Caelo	Lidocaine	15212002	3
Caelo	Lidocaine	15212001	1
Caelo	Lidocaine	13352415	3
Caelo	Lidocaine	15212001	1
Caelo	Lidocaine	160877708	1
Caelo	Lidocaine	16087708	1
Caelo	Lidocaine	16289504	3
Caelo	Lidocaine	17348503	7
Caelo	Lidocaine	12163811	2
Caelo	Lidocaine	15212003	3
Caelo	Lidocaine	14299305	2
Caelo	Lidocaine	13352405	1
Caelo	Lidocaine	17G03B08341308	1
Caelo	Lidocaine	12163813	4
Caelo	Lidocaine	13352411	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lidocaine	12294405	1
Caelo	Lidocaine	13352411	1
Caelo	Lidocaine	13J30-B02-292221	5
Caelo	Lidocaine	133524	1
Caelo	Lidocaine	16289502	9
Caelo	Lidocaine	13352413	3
Caelo	Lidocaine	14299301	1
Caelo	Lidocaine	14299303	2
Caelo	Lidocaine	17183103	2
Caelo	Lidocaine	13352406	1
Caelo	Lidocaine	13352403	1
Caelo	Lidocaine	122294405	1
Caelo	Lidocaine	13352414	2
Caelo	Lidocaine	9205919	1
Caelo	Lidocaine	12294405	1
Caelo	Lidocaine	12294408	2
Caelo	Lidocaine	12294412	9
Caelo	Lidocaine	13352408	8
Caelo	Lidocaine	12163802	2
Caelo	Lidocaine	16087711	8
Caelo	Lidocaine	16087704	7
Caelo	Lidocaine	15332314	4
Caelo	Lidocaine	16F21-806-333787	1
Caelo	Lidocaine	16087707	7
Caelo	Lidocaine	16087704	1
Caelo	Lidocaine	1804404	1
Caelo	Lidocaine	15332303	1
Caelo	Lidocaine	1629508	1
Caelo	Lidocaine	17183102	2
Caelo	Lidocaine	17183108	6
Caelo	Lidocaine	17183105	1
Caelo	Lidocaine	16087710	2
Caelo	Lidocaine	9205916	2
Caelo	Lidocaine	16087709	1
Caelo	Lidocaine	18044413	3
Caelo	Lidocaine	1804410	1
Caelo	Lidocaine	17K23-B04-353247	3
Caelo	Lidocaine	17183106	1
Caelo	Lidocaine	1612021-04	2
Caelo	Lidocaine	18044402	1
Caelo	Lidocaine	18044404	1
Caelo	Lidocaine	17183109	4
Caelo	Lidocaine	1718308	1
Caelo	Lidocaine	1783109	1
Caelo	Lidocaine	16087702	1
Caelo	Lidocaine	15332310	3
Caelo	Lidocaine	17348502	3
Caelo	Lidocaine	15332302	3
Caelo	Lidocaine	15H-B07-318310	1
Caelo	Lidocaine	16289508	13
Caelo	Lidocaine	12163811	1
Caelo	Lidocaine	16087712	6
Caelo	Lidocaine	14299301	1
Caelo	Lidocaine	1512003	1
Caelo	Lidocaine	15332306	3
Caelo	Lidocaine	5111201	1
Caelo	Lidocaine	16289507	2
Caelo	Lidocaine	12163804	2
Caelo	Lidocaine	12294412	2
Caelo	Lidocaine	16087703	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lidocaine	15332309	2
Caelo	Lidocaine	8629595	1
Caelo	Lidocaine	12163803	3
Caelo	Lidocaine hydrochlori...	17292304	13
Caelo	Lidocaine hydrochlori...	4111703	1
Caelo	Lidocaine hydrochlori...	17204803	7
Caelo	Lidocaine hydrochlori...	6121706	1
Caelo	Lidocaine hydrochlori...	812-16001	1
Caelo	Lidocaine hydrochlori...	3071711	1
Caelo	Lidocaine hydrochlori...	160602403	1
Caelo	Lidocaine hydrochlori...	16222208	6
Caelo	Lidocaine hydrochlori...	168020-0003	1
Caelo	Lidocaine hydrochlori...	16222214	15
Caelo	Lidocaine hydrochlori...	1606024	1
Caelo	Lidocaine hydrochlori...	152951602714117	1
Caelo	Lidocaine hydrochlori...	16222203	7
Caelo	Lidocaine hydrochlori...	1606024-02	49
Caelo	Lidocaine hydrochlori...	16222202	2
Caelo	Lidocaine hydrochlori...	15299516	6
Caelo	Lidocaine hydrochlori...	15299513	9
Caelo	Lidocaine hydrochlori...	15299514	6
Caelo	Lidocaine hydrochlori...	17204804	23
Caelo	Lidocaine hydrochlori...	1517-B01-315828	1
Caelo	Lidocaine hydrochlori...	15299509	2
Caelo	Lidocaine hydrochlori...	16083001	1
Caelo	Lidocaine hydrochlori...	15215312	8
Caelo	Lidocaine hydrochlori...	15299506	3
Caelo	Lidocaine hydrochlori...	160719LC	1
Caelo	Lidocaine hydrochlori...	15117-B01	1
Caelo	Lidocaine hydrochlori...	1230817/160602402	1
Caelo	Lidocaine hydrochlori...	13343011	3
Caelo	Lidocaine hydrochlori...	15299503	9
Caelo	Lidocaine hydrochlori...	26071303	1
Caelo	Lidocaine hydrochlori...	18179911	2
Caelo	Lidocaine hydrochlori...	2021903	1
Caelo	Lidocaine hydrochlori...	13118802	2
Caelo	Lidocaine hydrochlori...	1729313	1
Caelo	Lidocaine hydrochlori...	13118807	9
Caelo	Lidocaine hydrochlori...	AKD14	1
Caelo	Lidocaine hydrochlori...	4111302	1
Caelo	Lidocaine hydrochlori...	LDH/118024	1
Caelo	Lidocaine hydrochlori...	28111812	1
Caelo	Lidocaine hydrochlori...	13(1808031-01)	1
Caelo	Lidocaine hydrochlori...	17292311	3
Caelo	Lidocaine hydrochlori...	23011403	1
Caelo	Lidocaine hydrochlori...	17292313	3
Caelo	Lidocaine hydrochlori...	13118812	3
Caelo	Lidocaine hydrochlori...	13118811	2
Caelo	Lidocaine hydrochlori...	26041401	1
Caelo	Lidocaine hydrochlori...	1912M-01775	1
Caelo	Lidocaine hydrochlori...	13343001	2
Caelo	Lidocaine hydrochlori...	1515I-01775	1
Caelo	Lidocaine hydrochlori...	17292307	1
Caelo	Lidocaine hydrochlori...	17292308	6
Caelo	Lidocaine hydrochlori...	13118808	4
Caelo	Lidocaine hydrochlori...	13343001Mg:25g	1
Caelo	Lidocaine hydrochlori...	12135511	3
Caelo	Lidocaine hydrochlori...	10320719	1
Caelo	Lidocaine hydrochlori...	10320703	1
Caelo	Lidocaine hydrochlori...	11302906	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lidocaine hydrochlori...	10320707	2
Caelo	Lidocaine hydrochlori...	11091203	1
Caelo	Lidocaine hydrochlori...	10320704	1
Caelo	Lidocaine hydrochlori...	25091202	1
Caelo	Lidocaine hydrochlori...	120203501	1
Caelo	Lidocaine hydrochlori...	12176101	8
Caelo	Lidocaine hydrochlori...	12135503	2
Caelo	Lidocaine hydrochlori...	5211E-01775	1
Caelo	Lidocaine hydrochlori...	12135508	5
Caelo	Lidocaine hydrochlori...	11302904	1
Caelo	Lidocaine hydrochlori...	12176102	8
Caelo	Lidocaine hydrochlori...	10320717	1
Caelo	Lidocaine hydrochlori...	12176103	9
Caelo	Lidocaine hydrochlori...	17292316	1
Caelo	Lidocaine hydrochlori...	11302902	2
Caelo	Lidocaine hydrochlori...	11302905	1
Caelo	Lidocaine hydrochlori...	2051301	1
Caelo	Lidocaine hydrochlori...	12135509	4
Caelo	Lidocaine hydrochlori...	18B23-B03-355046	1
Caelo	Lidocaine hydrochlori...	12176109	2
Caelo	Lidocaine hydrochlori...	12135504	3
Caelo	Lidocaine hydrochlori...	13343004	7
Caelo	Lidocaine hydrochlori...	1792313	1
Caelo	Lidocaine hydrochlori...	280501	1
Caelo	Lidocaine hydrochlori...	15215311	6
Caelo	Lidocaine hydrochlori...	13343015	2
Caelo	Lidocaine hydrochlori...	15215310	6
Caelo	Lidocaine hydrochlori...	1814U-01775	1
Caelo	Lidocaine hydrochlori...	18041501	1
Caelo	Lidocaine hydrochlori...	13343003	5
Caelo	Lidocaine hydrochlori...	14254001	3
Caelo	Lidocaine hydrochlori...	13343019	4
Caelo	Lidocaine hydrochlori...	14110101	1
Caelo	Lidocaine hydrochlori...	13343017	6
Caelo	Lidocaine hydrochlori...	13343016	1
Caelo	Lidocaine hydrochlori...	29081503	1
Caelo	Lidocaine hydrochlori...	15215301	2
Caelo	Lidocaine hydrochlori...	15215310	3
Caelo	Lidocaine hydrochlori...	14J20-B05-303133	3
Caelo	Lidocaine hydrochlori...	15215304	4
Caelo	Lidocaine hydrochlori...	14J20-B05-303131	1
Caelo	Lidocaine hydrochlori...	1521304	1
Caelo	Lidocaine hydrochlori...	13343004	1
Caelo	Lidocaine hydrochlori...	15215301	5
Caelo	Lidocaine hydrochlori...	15215302	1
Caelo	Lidocaine hydrochlori...	15215302	1
Caelo	Lidocaine hydrochlori...	12021503	1
Caelo	Lidocaine hydrochlori...	26051802	1
Caelo	Lidocaine hydrochlori...	17292310	1
Caelo	Lidocaine hydrochlori...	17292302	4
Caelo	Lidocaine hydrochlori...	13343008	6
Caelo	Lidocaine hydrochlori...	24071401	1
Caelo	Lidocaine hydrochlori...	13343007	12
Caelo	Lidocaine hydrochlori...	15215309	3
Caelo	Lidocaine hydrochlori...	1403025-02	15
Caelo	Lidocaine hydrochlori...	15215308	3
Caelo	Lidocaine hydrochlori...	16101407	1
Caelo	Lidocaine hydrochlori...	13343012	16
Caelo	Lidocaine hydrochlori...	4051703	1
Caelo	Lidocaine hydrochlori...	17204810	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lidocaine hydrochlori...	13343014	5
Caelo	Lidocaine hydrochlori...	5113E-01775	1
Caelo	Lime flower, cut	13101701	7
Caelo	Lime flower, cut	12170702	2
Caelo	Lime flower, cut	11086307	3
Caelo	Lime flower, cut	13337403	5
Caelo	Lime flower, cut	13337403-348012102	1
Caelo	Lime flower, cut	13101708	1
Caelo	Lime flower, cut	15068302	1
Caelo	Lime flower, cut	15111301	1
Caelo	Lime flower, cut	Int.Ch:210140	1
Caelo	Lime flower, cut	272958	7
Caelo	Lime flower, cut	12170703	1
Caelo	Lime flower, cut	1210702	1
Caelo	Lime flower, cut	12170701	1
Caelo	Lime flower, cut	11086305	7
Caelo	Lime flower, cut	10273606	1
Caelo	Lime flower, cut	268709	2
Caelo	Lime flower, cut	9206604	1
Caelo	Lime flower, cut	28031302	1
Caelo	Lime flower, cut	13101702	7
Caelo	Lime flower, cut	12170702-3355319071	1
Caelo	Lime flower, cut	10273608	2
Caelo	Lime leaf, cut	1611A121120	1
Caelo	Linseed	11199006	1
Caelo	Linseed	16355402	3
Caelo	Linseed	12321601	1
Caelo	Linseed	13262704	2
Caelo	Linseed	12092104	2
Caelo	Linseed	12382002	2
Caelo	Linseed	28051303	1
Caelo	Linseed	12321602	4
Caelo	Linseed	12321604	2
Caelo	Linseed	1232006	1
Caelo	Linseed	12382004	1
Caelo	Linseed	12382006	3
Caelo	Linseed	12321607	1
Caelo	Linseed	18111302	1
Caelo	Linseed	13296705	5
Caelo	Linseed	13412804	1
Caelo	Linseed	13140305	1
Caelo	Linseed	12321606	2
Caelo	Linseed	7370A131104	1
Caelo	Linseed	1326703	2
Caelo	Linseed	28051403	1
Caelo	Linseed	13262703	3
Caelo	Linseed	13296702	3
Caelo	Linseed	12092103	3
Caelo	Linseed	13262702	1
Caelo	Linseed	11365507	1
Caelo	Linseed	7307A141201	2
Caelo	Linseed	17067204	1
Caelo	Linseed	7307A170619	2
Caelo	Linseed	21081703	1
Caelo	Linseed	16355401	2
Caelo	Linseed	16040005	1
Caelo	Linseed	16040002	1
Caelo	Linseed	16011403	3
Caelo	Linseed	16040001	1
Caelo	Linseed	3615A-02240	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Linseed	11365508	1
Caelo	Linseed	15321902	1
Caelo	Linseed	15006605	1
Caelo	Linseed	15006005	1
Caelo	Linseed	15061201	1
Caelo	Linseed	12092102	4
Caelo	Linseed	15006602	1
Caelo	Linseed	13296707	1
Caelo	Linseed	15053702	1
Caelo	Linseed	15006601	2
Caelo	Linseed	12186901	2
Caelo	Linseed	11199003	1
Caelo	Linseed	15006607	1
Caelo	Linseed	4713U02240	1
Caelo	Liquorice root, peeled	12227112	1
Caelo	Liquorice root, peeled	18121302	1
Caelo	Liquorice root, peeled	10011301	1
Caelo	Liquorice root, peeled	13194101	5
Caelo	Liquorice root, peeled	13104504	4
Caelo	Liquorice root, peeled	13104503	3
Caelo	Liquorice root, peeled	277004	1
Caelo	Liquorice root, peeled	638a120806	1
Caelo	Liquorice root, peeled	11351910	2
Caelo	Liquorice root, peeled	6381A111004	14
Caelo	Liquorice root, peeled	1135904	1
Caelo	Liquorice root, peeled	11351913	1
Caelo	Liquorice root, peeled	1351911	1
Caelo	Liquorice root, peeled	11351911	4
Caelo	Liquorice root, peeled	11351905	3
Caelo	Liquorice root, natur...	13037904	1
Caelo	Liquorice root, natur...	12267103	1
Caelo	Liquorice root, natur...	13037902	1
Caelo	Liquorice root, natur...	12267104	2
Caelo	Liquorice root, natur...		1
Caelo	Liquorice root, natur...	10110912	1
Caelo	Lovage root, cut	12169012	1
Caelo	Lovage root, cut	12169007	1
Caelo	Lovage root, cut	12169008	4
Caelo	Lovage root, cut	12169006	3
Caelo	Lovage root, cut	12169010	2
Caelo	Lovage root, cut	12169009	1
Caelo	Lovage root, cut	12169002	2
Caelo	Lovage root, cut	107063	1
Caelo	Lovage root, cut	11181611	3
Caelo	Lovage root, cut	4312M-02150	1
Caelo	Lovage root, cut	11181609	3
Caelo	Lovage root, cut	12189001	1
Caelo	Lovage root, cut	13385005	1
Caelo	Lovage root, cut	12189002	1
Caelo	L-tryptophan	131046	1
Caelo	Lungwort herb, cut	11021709	2
Caelo	Lungwort herb, cut	11021707	1
Caelo	Lungwort herb, cut	13029302	24
Caelo	Lungwort herb, cut	4331b111205	1
Caelo	Lungwort herb, cut	13355801	2
Caelo	Lungwort herb, cut	13029306	3
Caelo	Lungwort herb, cut	11134201	1
Caelo	Lungwort herb, cut	11134205	3
Caelo	Lungwort herb, cut	11134203	12
Caelo	Lungwort herb, cut	11134202	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Lycopodium herb, cut	11255005	3
Caelo	Lycopodium herb, cut	11255008	12
Caelo	Lycopodium herb, cut	13101501	1
Caelo	Lycopodium herb, cut	1310502	1
Caelo	Lycopodium herb, cut	11255018	1
Caelo	Lycopodium herb, cut	11255010	1
Caelo	Lycopodium herb, cut	11255017	6
Caelo	Lycopodium herb, cut	11255012	2
Caelo	Lycopodium herb, cut	11255016	2
Caelo	Lycopodium herb, cut	11255014	7
Caelo	Lycopodium herb, cut	11255011	5
Caelo	Lycopodium herb, cut	11255003	1
Caelo	Lycopodium herb, cut	11255009	2
Caelo	Lycopodium herb, cut	9378812	1
Caelo	Macrogols 1500	18129505	1
Caelo	Macrogols 1500	18045703	1
Caelo	Macrogols 1500	18129506	1
Caelo	Macrogols 1500	14072412	1
Caelo	Macrogols 1500	17120001	1
Caelo	Macrogols 1500	13140906	1
Caelo	Macrogols 1500	17120002	2
Caelo	Macrogols 1500	10105901	1
Caelo	Macrogols 1500	12109409	1
Caelo	Macrogols 1500	15413111	1
Caelo	Macrogols 1500	13140901	1
Caelo	Macrogols 1500	13140902	2
Caelo	Macrogols 1500	16138407	1
Caelo	Macrogols 1500	15413105	5
Caelo	Macrogols 1500	15413101	1
Caelo	Macrogols 1500	15120007	1
Caelo	Macrogols 1500	18129503	1
Caelo	Macrogols 1500	11071309	1
Caelo	Macrogols 4000	13136010	1
Caelo	Macrogols 4000	18106109	2
Caelo	Macrogols 4000	18106106	1
Caelo	Macrogols 4000	18106103	4
Caelo	Macrogols 4000	15225013	1
Caelo	Macrogols 4000	18106104	3
Caelo	Macrogols 4000	9061805	1
Caelo	Macrogols 4000	15376503	3
Caelo	Macrogols 4000	16023903	4
Caelo	Macrogols 4000	18106101	2
Caelo	Macrogols 4000	16152001	13
Caelo	Macrogols 4000	17337902/1300518	1
Caelo	Macrogols 4000	16152003	4
Caelo	Macrogols 4000	15376515	1
Caelo	Macrogols 4000	17337902/2190218	1
Caelo	Macrogols 4000	17337902/1250418	1
Caelo	Macrogols 4000	30041803	1
Caelo	Macrogols 4000	17303901	7
Caelo	Macrogols 4000	160912M	1
Caelo	Macrogols 4000	161017M/15406304	1
Caelo	Macrogols 4000	161202M/16152004	1
Caelo	Macrogols 4000	16023904	5
Caelo	Macrogols 4000	15376518	1
Caelo	Macrogols 4000	4123585	1
Caelo	Macrogols 4000	15225001	1
Caelo	Macrogols 4000	18215403	3
Caelo	Macrogols 4000	18244803	2
Caelo	Macrogols 4000	18244805	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 4000	18244802	2
Caelo	Macrogols 4000	18215404	2
Caelo	Macrogols 4000	17337901	5
Caelo	Macrogols 4000	16308904	7
Caelo	Macrogols 4000	190920	1
Caelo	Macrogols 4000	15159515	5
Caelo	Macrogols 4000	15159514	4
Caelo	Macrogols 4000	15335001	9
Caelo	Macrogols 4000	15159516	1
Caelo	Macrogols 4000	15335003	6
Caelo	Macrogols 4000	15335007	4
Caelo	Macrogols 4000	18106105	10
Caelo	Macrogols 4000	15225003	4
Caelo	Macrogols 4000	15406301	9
Caelo	Macrogols 4000	15335002	3
Caelo	Macrogols 4000	18106110	2
Caelo	Macrogols 4000	15406303	3
Caelo	Macrogols 4000	17337906	1
Caelo	Macrogols 4000	39393137	1
Caelo	Macrogols 4000	16267305	8
Caelo	Macrogols 4000	17035606	5
Caelo	Macrogols 4000	17119001	2
Caelo	Macrogols 4000	16358106	2
Caelo	Macrogols 4000	17035608	1
Caelo	Macrogols 4000	170817M/17035603	1
Caelo	Macrogols 4000	16358101	2
Caelo	Macrogols 4000	170710M/17035603	1
Caelo	Macrogols 4000	1200617/16368105	1
Caelo	Macrogols 4000	17035604	3
Caelo	Macrogols 4000	17035603	1
Caelo	Macrogols 4000	170612M/17035603	1
Caelo	Macrogols 4000	16308907	13
Caelo	Macrogols 4000	170529M/16267308	1
Caelo	Macrogols 4000	170504M/16308906	1
Caelo	Macrogols 4000	1030517	1
Caelo	Macrogols 4000	5041701	1
Caelo	Macrogols 4000	170324M/16267308	1
Caelo	Macrogols 4000	16308905	3
Caelo	Macrogols 4000	16267308/170227M	1
Caelo	Macrogols 4000	170928M/17119002	1
Caelo	Macrogols 4000	16267311	6
Caelo	Macrogols 4000	16308908	3
Caelo	Macrogols 4000	17119005	2
Caelo	Macrogols 4000	170105M16152004	1
Caelo	Macrogols 4000	17303902	9
Caelo	Macrogols 4000	17208106	1
Caelo	Macrogols 4000	15376519	1
Caelo	Macrogols 4000	17303903	1
Caelo	Macrogols 4000	16267307	11
Caelo	Macrogols 4000	1720103	1
Caelo	Macrogols 4000	143526	1
Caelo	Macrogols 4000	16308902	9
Caelo	Macrogols 4000	16358107	1
Caelo	Macrogols 4000	17208105	1
Caelo	Macrogols 4000	16358107/1111217	1
Caelo	Macrogols 4000	171211M/17119003	1
Caelo	Macrogols 4000	17208104	7
Caelo	Macrogols 4000	28111706	1
Caelo	Macrogols 4000	1720813	1
Caelo	Macrogols 4000	1131117/16358107	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 4000	1260917/16358105	1
Caelo	Macrogols 4000	171103M/17119002	1
Caelo	Macrogols 4000	17208103	13
Caelo	Macrogols 4000	14J03-B07-304700	3
Caelo	Macrogols 4000	16267308	3
Caelo	Macrogols 4000	13335102	1
Caelo	Macrogols 4000	13335107	1
Caelo	Macrogols 4000	13356621	1
Caelo	Macrogols 4000	18315304	2
Caelo	Macrogols 4000	13356616	2
Caelo	Macrogols 4000	14061705	5
Caelo	Macrogols 4000	123Test	2
Caelo	Macrogols 4000		1
Caelo	Macrogols 4000	12074114	1
Caelo	Macrogols 4000	13335120	1
Caelo	Macrogols 4000	14180503	1
Caelo	Macrogols 4000	13356601	2
Caelo	Macrogols 4000	14061708	4
Caelo	Macrogols 4000	1180513	1
Caelo	Macrogols 4000	1401708	1
Caelo	Macrogols 4000	1234	1
Caelo	Macrogols 4000	14180511	2
Caelo	Macrogols 4000	14061708	5
Caelo	Macrogols 4000	13356617	6
Caelo	Macrogols 4000	13356611	2
Caelo	Macrogols 4000	13356602	3
Caelo	Macrogols 4000	1316003	1
Caelo	Macrogols 4000	12361518	3
Caelo	Macrogols 4000	IntCh217002	1
Caelo	Macrogols 4000	3940608	1
Caelo	Macrogols 4000	12318019	1
Caelo	Macrogols 4000	13136001	3
Caelo	Macrogols 4000	12318011	1
Caelo	Macrogols 4000	12361513	3
Caelo	Macrogols 4000	12361506	3
Caelo	Macrogols 4000	1010023-02	1
Caelo	Macrogols 4000	12074102	1
Caelo	Macrogols 4000	12149113	2
Caelo	Macrogols 4000	13136011	4
Caelo	Macrogols 4000	12149108	1
Caelo	Macrogols 4000	13335106	1
Caelo	Macrogols 4000	13356601	2
Caelo	Macrogols 4000	11172303	1
Caelo	Macrogols 4000	14196007	1
Caelo	Macrogols 4000	13136009	3
Caelo	Macrogols 4000	14061713	1
Caelo	Macrogols 4000	14180524	1
Caelo	Macrogols 4000	18244809	1
Caelo	Macrogols 4000	18315301	5
Caelo	Macrogols 4000	15159502	3
Caelo	Macrogols 4000	18315303	2
Caelo	Macrogols 4000	14196017	2
Caelo	Macrogols 4000	13081501	1
Caelo	Macrogols 4000	18315305	1
Caelo	Macrogols 4000	15159504	2
Caelo	Macrogols 4000	14352613	1
Caelo	Macrogols 4000	14180520	4
Caelo	Macrogols 4000	12074109	3
Caelo	Macrogols 4000	1172404	1
Caelo	Macrogols 4000	23061505	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Macrogols 4000	12318027	1
Caelo	Macrogols 4000	14352606	2
Caelo	Macrogols 4000	14180504	1
Caelo	Macrogols 4000	14352614	1
Caelo	Macrogols 4000	14352613	4
Caelo	Macrogols 4000	12074116	3
Caelo	Macrogols 4000	14352611	2
Caelo	Macrogols 4000	14180508	1
Caelo	Macrogols 4000	14180525	3
Caelo	Macrogols 4000	14180522	1
Caelo	Macrogols 4000	14180508	1
Caelo	Macrogols 4000	14196013	1
Caelo	Macrogols 4000	14352608	1
Caelo	Macrogols 4000	14352618	1
Caelo	Macrogols 4000	14352607	2
Caelo	Magnesium carbonate, ...	7011404	1
Caelo	Magnesium carbonate, ...	2091303	1
Caelo	Magnesium carbonate, ...	13123505	1
Caelo	Magnesium carbonate, ...	13123502	1
Caelo	Magnesium carbonate, ...	12135404	1
Caelo	Magnesium carbonate, ...	13123504	2
Caelo	Magnesium chloride he...	12354615	8
Caelo	Magnesium chloride he...	12354606	4
Caelo	Magnesium chloride he...	13194706	3
Caelo	Magnesium chloride he...	12354610	2
Caelo	Magnesium chloride he...	1512I-01828	1
Caelo	Magnesium chloride he...	13194702	3
Caelo	Magnesium chloride he...	12141807	3
Caelo	Magnesium chloride he...	12354618	2
Caelo	Magnesium chloride he...	12141806	2
Caelo	Magnesium chloride he...	12354613	4
Caelo	Magnesium chloride he...	12354608	1
Caelo	Magnesium chloride he...	12141811	2
Caelo	Magnesium chloride he...	12141817	2
Caelo	Magnesium chloride he...	12354605	4
Caelo	Magnesium chloride he...	12354612	7
Caelo	Magnesium chloride he...	11345623	2
Caelo	Magnesium chloride he...	11345606	1
Caelo	Magnesium chloride he...	1108024-01	1
Caelo	Magnesium chloride he...	11345609	1
Caelo	Magnesium chloride he...	12141814	1
Caelo	Magnesium chloride he...	5210I-01828	1
Caelo	Magnesium chloride he...	12354604	4
Caelo	Magnesium chloride he...	12141808	1
Caelo	Magnesium chloride he...	11345624	1
Caelo	Magnesium chloride he...	17034304	1
Caelo	Magnesium chloride he...	16094014	1
Caelo	Magnesium chloride he...	16094004	1
Caelo	Magnesium chloride he...	13194721	2
Caelo	Magnesium chloride he...	1304C004	2
Caelo	Magnesium chloride he...	1.02.14	1
Caelo	Magnesium chloride he...	13194707	4
Caelo	Magnesium chloride he...	13194713	1
Caelo	Magnesium chloride he...	13194713-3490312103	1
Caelo	Magnesium chloride he...	13194701	1
Caelo	Magnesium chloride he...	12141803	1
Caelo	Magnesium chloride he...	12354609	7
Caelo	Magnesium hydrogen ci...	16E02-B03-321485	1
Caelo	Magnesium hydrogen ci...	13233913	6
Caelo	Magnesium hydrogen ci...	90396169	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Magnesium hydrogen ci...	13233905	1
Caelo	Magnesium hydrogen ci...	13233911	3
Caelo	Magnesium hydrogen ci...	411052-2-P	1
Caelo	Magnesium hydrogen ci...	133233907	1
Caelo	Magnesium hydrogen ci...	11021401	1
Caelo	Magnesium hydrogen ci...	913M-01829	1
Caelo	Magnesium hydrogen ci...	13389103	8
Caelo	Magnesium hydrogen ci...	13389106	2
Caelo	Magnesium hydrogen ci...	11061403	1
Caelo	Magnesium hydrogen ci...	17061401	1
Caelo	Magnesium hydrogen ci...	13233910	5
Caelo	Magnesium hydrogen ci...	27111306	1
Caelo	Magnesium hydrogen ci...	13389111	7
Caelo	Magnesium hydrogen ci...	13389107	2
Caelo	Magnesium hydrogen ci...	14I12-B02-298854	2
Caelo	Magnesium hydrogen ci...	14265503	6
Caelo	Magnesium hydrogen ci...	7011504	1
Caelo	Magnesium hydrogen ci...	1014E-01829	1
Caelo	Magnesium hydrogen ci...	14265505	2
Caelo	Magnesium hydrogen ci...	14265511	1
Caelo	Magnesium hydrogen ci...	14265504	5
Caelo	Magnesium hydrogen ci...	14265514Gehe14.04.	1
Caelo	Magnesium hydrogen ci...	14265514	1
Caelo	Magnesium hydrogen ci...	14265517	1
Caelo	Magnesium hydrogen ci...	14265513	2
Caelo	Magnesium hydrogen ci...	133389111	1
Caelo	Magnesium hydrogen ci...	14265521	2
Caelo	Magnesium hydrogen ci...	13233904	5
Caelo	Magnesium hydrogen ci...	13233901	9
Caelo	Magnesium hydrogen ci...	11372902	1
Caelo	Magnesium hydrogen ci...	12175801	1
Caelo	Magnesium hydrogen ci...	12175807	3
Caelo	Magnesium hydrogen ci...	12175802	2
Caelo	Magnesium hydrogen ci...	26101203	1
Caelo	Magnesium hydrogen ci...	12175810	5
Caelo	Magnesium hydrogen ci...	12238002	2
Caelo	Magnesium hydrogen ci...	12238012	1
Caelo	Magnesium hydrogen ci...	11044103	2
Caelo	Magnesium hydrogen ci...	12238004	3
Caelo	Magnesium hydrogen ci...	12368003	9
Caelo	Magnesium hydrogen ci...	121033-01	1
Caelo	Magnesium hydrogen ci...	11372901	1
Caelo	Magnesium hydrogen ci...	12368005	1
Caelo	Magnesium hydrogen ci...	11193104	1
Caelo	Magnesium hydrogen ci...	19061316	1
Caelo	Magnesium hydrogen ci...	12368011	8
Caelo	Magnesium hydrogen ci...	80446188	1
Caelo	Magnesium hydrogen ci...	2368011	1
Caelo	Magnesium hydrogen ci...	12368002	3
Caelo	Magnesium hydrogen ci...	6081301	1
Caelo	Magnesium hydrogen ci...	12368017	2
Caelo	Magnesium hydrogen ci...	12091303	1
Caelo	Magnesium hydrogen ci...	13194601	4
Caelo	Magnesium hydrogen ci...	13233906	9
Caelo	Magnesium hydrogen ci...	13194603	1
Caelo	Magnesium hydrogen ci...	12368004	2
Caelo	Magnesium hydrogen ci...	1214	1
Caelo	Magnesium hydrogen ci...	14265516	3
Caelo	Magnesium hydrogen ci...	15296801	3
Caelo	Magnesium hydrogen ci...	15296702	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Magnesium hydrogen ci...	18216803	1
Caelo	Magnesium hydrogen ci...	17353104	2
Caelo	Magnesium hydrogen ci...	17173406	2
Caelo	Magnesium hydrogen ci...	17173404	1
Caelo	Magnesium hydrogen ci...	17173401	3
Caelo	Magnesium hydrogen ci...	27101704	1
Caelo	Magnesium hydrogen ci...	16341710	8
Caelo	Magnesium hydrogen ci...	16341711	2
Caelo	Magnesium hydrogen ci...	16341709	4
Caelo	Magnesium hydrogen ci...	16341707	1
Caelo	Magnesium hydrogen ci...	15415014	5
Caelo	Magnesium hydrogen ci...	15415013	2
Caelo	Magnesium hydrogen ci...	1415013	1
Caelo	Magnesium hydrogen ci...	16341708	3
Caelo	Magnesium hydrogen ci...	15415009	5
Caelo	Magnesium hydrogen ci...	15415010	10
Caelo	Magnesium hydrogen ci...	315E-01829	1
Caelo	Magnesium hydrogen ci...	1601005-01	15
Caelo	Magnesium hydrogen ci...	15296712	2
Caelo	Magnesium hydrogen ci...	15296709	4
Caelo	Magnesium hydrogen ci...	15415001	5
Caelo	Magnesium hydrogen ci...	15415003	6
Caelo	Magnesium hydrogen ci...	15296701	8
Caelo	Magnesium hydrogen ci...	15415002	1
Caelo	Magnesium hydrogen ci...	15415005	4
Caelo	Magnesium hydrogen ci...	14265520	1
Caelo	Magnesium hydrogen ci...	516Q-01829	1
Caelo	Magnesium hydrogen ci...	15415004	2
Caelo	Magnesium orotate dihydro...	14093607	2
Caelo	Magnesium stearate	15305307	1
Caelo	Magnesium stearate	15305302	1
Caelo	Magnesium stearate	15305323	1
Caelo	Magnesium stearate	15305315	3
Caelo	Magnesium stearate	13095515	1
Caelo	Magnesium stearate	16206602	2
Caelo	Magnesium stearate	13095536	1
Caelo	Magnesium stearate	12040307	1
Caelo	Magnesium stearate	15155001	1
Caelo	Magnesium stearate	16206610	1
Caelo	Magnesium sulfate hep...	6012014B	1
Caelo	Magnesium sulfate hep...	1213E-01835	1
Caelo	Magnesium sulfate hep...	11369252	1
Caelo	Magnesium sulfate hep...	12c27-n06	1
Caelo	Magnesium sulfate hep...	13183702	2
Caelo	Magnesium sulfate hep...	12119901/2	1
Caelo	Magnesium sulfate hep...	12336915	3
Caelo	Magnesium sulfate hep...	12336926	1
Caelo	Magnesium sulfate hep...	13030104	1
Caelo	Magnesium sulfate hep...	1206019-01	5
Caelo	Magnesium sulfate hep...	10104201	1
Caelo	Magnesium sulfate hep...	10371111	1
Caelo	Magnesium sulfate hep...	11369526	13
Caelo	Magnesium sulfate hep...	11369529	1
Caelo	Magnesium sulfate hep...	11194001	1
Caelo	Magnesium sulfate hep...	11369525	6
Caelo	Magnesium sulfate hep...	11369509	1
Caelo	Magnesium sulfate hep...	90744129	1
Caelo	Magnesium sulfate hep...	12336914	5
Caelo	Magnesium sulfate hep...	12138521	4
Caelo	Magnesium sulfate hep...	13030107	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Magnesium sulfate hep...	113695625	1
Caelo	Magnesium sulfate, dr...	11265902	1
Caelo	Magnesium sulfate, dr...	11265908	2
Caelo	Magnesium sulfate, dr...	10213213	2
Caelo	Magnesium sulfate, dr...	11265903	1
Caelo	Magnesium sulfate, dr...	12157202	1
Caelo	Magnesium sulfate, dr...	M-220413-2	1
Caelo	Magnesium sulfate, dr...	12157207	1
Caelo	Magnesium sulfate, dr...	12157208	1
Caelo	Maize style, cut	13277201	1
Caelo	Maize style, cut	13277207	1
Caelo	Maize style, cut	11280404	2
Caelo	Maize style, cut	13277202	1
Caelo	Maize style, cut	10325902	2
Caelo	Maize style, cut	10325906	7
Caelo	Maize style, cut	10325904	1
Caelo	Maize style, cut	10325907	4
Caelo	Maize style, cut	10325905	1
Caelo	Maize style, cut	11280405	10
Caelo	Maize style, cut	280363	1
Caelo	Maize style, cut	13277205	1
Caelo	Male speedwell wort, ...	3051306	1
Caelo	Male speedwell wort, ...	28061303	1
Caelo	Male speedwell wort, ...	12162505	2
Caelo	Male speedwell wort, ...	3091305	1
Caelo	Male speedwell wort, ...	13166701	3
Caelo	Male speedwell wort, ...	14111303	1
Caelo	Male speedwell wort, ...	11265503	5
Caelo	Male speedwell wort, ...	4661c111011	1
Caelo	Male speedwell wort, ...	1216504	1
Caelo	Male speedwell wort, ...	11265502	1
Caelo	Male speedwell wort, ...	12162503	1
Caelo	Male speedwell wort, ...	12162504	12
Caelo	Male speedwell wort, ...	20111303	1
Caelo	Male speedwell wort, ...	11061305	1
Caelo	Mallow flower, cut	13351307	1
Caelo	Mallow flower, cut	13351309	1
Caelo	Mallow flower, cut	18111303	1
Caelo	Mallow flower, cut	13267303	2
Caelo	Mallow flower, cut	12187501-3414310132	1
Caelo	Mallow flower, cut	13267301-3409310112	1
Caelo	Mallow flower, cut	13040402	1
Caelo	Mallow flower, cut	9311104	1
Caelo	Mallow flower, cut	10357902	1
Caelo	Mallow flower, cut	13267302	1
Caelo	Mallow flower, cut	20061201	1
Caelo	Mallow flower, cut	8121202	1
Caelo	Mallow flower, cut	27021301	1
Caelo	Mallow flower, cut	12187405	4
Caelo	Mallow flower, cut	12187404	4
Caelo	Mallow flower, cut	12187501	6
Caelo	Mallow flower, cut	10357904	1
Caelo	Mallow flower, cut	11040106	4
Caelo	Mallow flower, cut	12187505	3
Caelo	Mallow flower, cut	13351301	2
Caelo	Mallow flower, cut	12187502	1
Caelo	Mallow flower, cut	13351303	2
Caelo	Mallow leaf, cut	28061203	1
Caelo	Mallow leaf, cut	10290515	1
Caelo	Mallow leaf, cut	12153403	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Mallow leaf, cut	2131B120529	3
Caelo	Mallow leaf, cut	12310001	2
Caelo	Mallow leaf, cut	12310004	4
Caelo	Mallow leaf, cut	10290502	1
Caelo	Mallow leaf, cut	12310006	1
Caelo	Mallow leaf, cut	12310008	1
Caelo	Mallow leaf, cut	20111302	1
Caelo	Mallow leaf, cut	13174207	1
Caelo	Mallow leaf, cut	13174213	1
Caelo	Mallow leaf, cut	14073902	1
Caelo	Mallow leaf, cut	12310009	2
Caelo	Mallow leaf, cut	13174217	1
Caelo	Manna, cut	15111738	1
Caelo	Manna, cut	15219010	3
Caelo	Manna, cut	11308901	2
Caelo	Manna, cut	11308909	2
Caelo	Manna, cut	11308907	2
Caelo	Manna, cut	11308911	2
Caelo	Manna, cut	25111301	1
Caelo	Manna, cut	13127003	5
Caelo	Manna, cut	11061404	1
Caelo	Manna, cut	25071403	1
Caelo	Manna, cut	11041401	1
Caelo	Manna, cut	14001303	3
Caelo	Manna, cut	14001310	3
Caelo	Manna, cut	15219003	2
Caelo	Manna, cut	15219005	1
Caelo	Manna, cut	13041712	1
Caelo	Manna, cut	15219008	2
Caelo	Manna, cut	28101404	1
Caelo	Manna, powdered	11003611	1
Caelo	Manna, powdered	11071311	1
Caelo	Manna, powdered	11308807	1
Caelo	Manna, powdered	14001402	1
Caelo	Manna, powdered	11308802	2
Caelo	Manna, powdered	14001408	1
Caelo	Manna, powdered	15101306	1
Caelo	Manna, powdered	15101312	1
Caelo	Mannitol	13189228	1
Caelo	Mannitol	14051208	2
Caelo	Mannitol	4414E-01844	1
Caelo	Mannitol	1400434-01	1
Caelo	Mannitol	13K20-B01-287704	1
Caelo	Mannitol	14051209	7
Caelo	Mannitol	1405127	1
Caelo	Mannitol	2101401	1
Caelo	Mannitol	14051207	1
Caelo	Mannitol	12305736	4
Caelo	Mannitol	14051217	2
Caelo	Mannitol	14051215	6
Caelo	Mannitol	14051205	6
Caelo	Mannitol	14051230	2
Caelo	Mannitol	30121417	1
Caelo	Mannitol	1423120	1
Caelo	Mannitol	14323126	4
Caelo	Mannitol	15184101	1
Caelo	Mannitol	14323115	4
Caelo	Mannitol	14323121	7
Caelo	Mannitol	14323114	3
Caelo	Mannitol	14323105	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Mannitol	14323103	2
Caelo	Mannitol	30031502	1
Caelo	Mannitol	14051229	5
Caelo	Mannitol	14323101	3
Caelo	Mannitol	13095013	1
Caelo	Mannitol	14051216	3
Caelo	Mannitol	14051223	1
Caelo	Mannitol	10011501	1
Caelo	Mannitol	14051224	2
Caelo	Mannitol	13189227	12
Caelo	Mannitol	13189208	6
Caelo	Mannitol	7071404	1
Caelo	Mannitol	9081303	1
Caelo	Mannitol	11003611	1
Caelo	Mannitol	12305709	16
Caelo	Mannitol	12305710	9
Caelo	Mannitol	5031304	1
Caelo	Mannitol	10100026	2
Caelo	Mannitol	9204719	1
Caelo	Mannitol	11219121	1
Caelo	Mannitol	12105438	2
Caelo	Mannitol	20151012-2	1
Caelo	Mannitol	11D01-N10	1
Caelo	Mannitol	12105421	9
Caelo	Mannitol	12105403	3
Caelo	Mannitol	11219124	2
Caelo	Mannitol	12105420	7
Caelo	Mannitol	12105405	1
Caelo	Mannitol	1110009-02	5
Caelo	Mannitol	12305725	12
Caelo	Mannitol	23082013C	1
Caelo	Mannitol	12305724	5
Caelo	Mannitol	12305726	1
Caelo	Mannitol	13189222	3
Caelo	Mannitol	16061401	1
Caelo	Mannitol	13189226	4
Caelo	Mannitol	13189207	8
Caelo	Mannitol	13189215	18
Caelo	Mannitol	14323127	3
Caelo	Mannitol	13189203-3566412081	1
Caelo	Mannitol	12305728	1
Caelo	Mannitol	14051204	8
Caelo	Mannitol	16011402	1
Caelo	Mannitol	13189203	8
Caelo	Mannitol	1230578	1
Caelo	Mannitol	12305740	2
Caelo	Mannitol	11308807	2
Caelo	Mannitol	12305725-3444311112	1
Caelo	Mannitol	13189201	8
Caelo	Mannitol	24101306	1
Caelo	Mannitol	12305734	2
Caelo	Mannitol	13189206	1
Caelo	Mannitol	18265411	1
Caelo	Mannitol	14323138	26
Caelo	Mannitol	14323125	2
Caelo	Mannitol	17I13-B02-340966	1
Caelo	Mannitol	17045009	7
Caelo	Mannitol	17045006	3
Caelo	Mannitol	17087909	1
Caelo	Mannitol	17087908	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Mannitol	18121701	1
Caelo	Mannitol	151117	1
Caelo	Mannitol	17045002	12
Caelo	Mannitol	16236603	1
Caelo	Mannitol	170450	1
Caelo	Mannitol	16200924	9
Caelo	Mannitol	17045004	2
Caelo	Mannitol	163403	1
Caelo	Mannitol	1610018-01	1
Caelo	Mannitol	16200912	3
Caelo	Mannitol	1071701	1
Caelo	Mannitol	17045003	8
Caelo	Mannitol	17045001	3
Caelo	Mannitol	17A12-B03	1
Caelo	Mannitol	17087916	14
Caelo	Mannitol	16L15-B04-329317	3
Caelo	Mannitol	200218P302	1
Caelo	Mannitol	17045008	5
Caelo	Mannitol	17258819	3
Caelo	Mannitol	18265502	4
Caelo	Mannitol	1809008-01	4
Caelo	Mannitol	18D09-B12-351437	1
Caelo	Mannitol	18D09-B13-352214	4
Caelo	Mannitol	17360112	5
Caelo	Mannitol	27101802	1
Caelo	Mannitol	11(16200912)	1
Caelo	Mannitol	17258812	3
Caelo	Mannitol	17360108	8
Caelo	Mannitol	17258808	3
Caelo	Mannitol	17258809	5
Caelo	Mannitol	17087915	2
Caelo	Mannitol	E911K	1
Caelo	Mannitol	9061808	1
Caelo	Mannitol	17J30-B10-3434414	1
Caelo	Mannitol	17360101	5
Caelo	Mannitol	170800501	1
Caelo	Mannitol	17087907	1
Caelo	Mannitol	17045010	3
Caelo	Mannitol	14323120	1
Caelo	Mannitol	16200916	2
Caelo	Mannitol	16200908	2
Caelo	Mannitol	1612101	1
Caelo	Mannitol	16212103	2
Caelo	Mannitol	15D28-B02-309192	4
Caelo	Mannitol	15197327	1
Caelo	Mannitol	15197312	1
Caelo	Mannitol	16D20-B03	2
Caelo	Mannitol	1519713	1
Caelo	Mannitol	15525	1
Caelo	Mannitol	15K25-B04	1
Caelo	Mannitol	48558	1
Caelo	Mannitol	2274635	1
Caelo	Mannitol	15197313	18
Caelo	Mannitol	1508012-02	14
Caelo	Mannitol	15197302	5
Caelo	Mannitol	15197301	13
Caelo	Mannitol	15197303	1
Caelo	Mannitol	14323124	4
Caelo	Mannitol	18265504	1
Caelo	Mannitol	14323129	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Mannitol	16212102	11
Caelo	Mannitol	16200913	3
Caelo	Mannitol	16014016	1
Caelo	Mannitol	160400801	1
Caelo	Mannitol	16200903	9
Caelo	Mannitol	16212106	4
Caelo	Mannitol	16340302	6
Caelo	Mannitol	8031704	1
Caelo	Mannitol	16212109	4
Caelo	Mannitol	16340301	14
Caelo	Mannitol	15197319	8
Caelo	Mannitol	16200901	2
Caelo	Mannitol	16200909	1
Caelo	Mannitol	16014005	7
Caelo	Mannitol	15197329	1
Caelo	Mannitol	16H29-B05-327813	2
Caelo	Mannitol	16200904	2
Caelo	Mannitol	1621212	1
Caelo	Mannitol	15197219	1
Caelo	Mannitol	16212104	2
Caelo	Mannitol	16G27-B03	1
Caelo	Mannitol	18265507	1
Caelo	Mannitol	162121	1
Caelo	Mannitol	16014014	1
Caelo	Mannitol	4516E-01844	1
Caelo	Mannitol	16014013	1
Caelo	Marian thistle fruit	27111203	1
Caelo	Marian thistle fruit	15409204	2
Caelo	Marian thistle fruit	13326202	3
Caelo	Marian thistle fruit	12189116	2
Caelo	Marian thistle fruit	18111301	1
Caelo	Marian thistle fruit	12189115	2
Caelo	Marian thistle fruit	12189111	3
Caelo	Marian thistle fruit	12189112	5
Caelo	Marian thistle fruit	12198107	1
Caelo	Marian thistle fruit	12189113	6
Caelo	Marian thistle fruit	12189109	2
Caelo	Marian thistle fruit	12189108	4
Caelo	Marian thistle fruit	12189105	1
Caelo	Marian thistle fruit	12189107	2
Caelo	Marian thistle fruit	12189104	2
Caelo	Marian thistle fruit	12189106	4
Caelo	Marian thistle fruit	1289112	1
Caelo	Marian thistle fruit	1126703	1
Caelo	Marian thistle fruit	11126705	1
Caelo	Marian thistle fruit	10270613	1
Caelo	Marian thistle fruit	10270614	2
Caelo	Marian thistle fruit	11126702	9
Caelo	Marian thistle fruit	1126705	1
Caelo	Marian thistle fruit	11126703	1
Caelo	Marian thistle fruit	11126701	1
Caelo	Marian thistle herb, ...	10315907	1
Caelo	Marian thistle herb, ...	17121204	1
Caelo	Marian thistle herb, ...	30101306	1
Caelo	Marian thistle herb, ...	10365605	1
Caelo	Marian thistle herb, ...	13130703	2
Caelo	Marian thistle herb, ...	intch:318080	1
Caelo	Marian thistle herb, ...	18041303	1
Caelo	Marian thistle herb, ...	10365601	4
Caelo	Marian thistle herb, ...	13130704	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Marian thistle herb, ...	62607127	1
Caelo	Marian thistle herb, ...	10365604	4
Caelo	Marian thistle herb, ...	13130701-3492411030	1
Caelo	Marian thistle herb, ...	10315909	1
Caelo	Marian thistle herb, ...	13386703-3543411013	1
Caelo	Marian thistle herb, ...	10315911	3
Caelo	Marian thistle herb, ...	13130701	2
Caelo	Marian thistle herb, ...	13130702	3
Caelo	Marian thistle herb, ...	1035911	1
Caelo	Marian thistle herb, ...	1313702	1
Caelo	Marjoram herb, rubbed	12168903	10
Caelo	Marjoram herb, rubbed	12168906	1
Caelo	Marjoram herb, rubbed	12168905	2
Caelo	Marjoram herb, rubbed	27011402	1
Caelo	Marjoram herb, rubbed	12168908	6
Caelo	Marjoram herb, rubbed	12168902	1
Caelo	Marjoram herb, rubbed	13113303	3
Caelo	Marjoram herb, rubbed	12168912	1
Caelo	Marjoram herb, rubbed	9292602	1
Caelo	Marjoram herb, rubbed	12168909	6
Caelo	Marjoram herb, rubbed	12168910	8
Caelo	Marjoram herb, rubbed	12168907	3
Caelo	Marjoram herb, rubbed	11212711	3
Caelo	Marjoram herb, rubbed	12168904	1
Caelo	Marjoram herb, rubbed	13113304	3
Caelo	Marjoram herb, rubbed	12168901	3
Caelo	Marjoram herb, rubbed	6051302	1
Caelo	Marrubium vulgare her...	12231604	4
Caelo	Marrubium vulgare her...	11253803	1
Caelo	Marrubium vulgare her...	313091	1
Caelo	Marrubium vulgare her...	12231603	1
Caelo	Marrubium vulgare her...	12231601	4
Caelo	Marrubium vulgare her...	11253802	1
Caelo	Marrubium vulgare her...	13289109	1
Caelo	Marrubium vulgare her...	1153803	1
Caelo	Marrubium vulgare her...	12231605	4
Caelo	Marshmallow leaf, cut	12328202	3
Caelo	Marshmallow leaf, cut	13199604	4
Caelo	Marshmallow leaf, cut	12328207	2
Caelo	Marshmallow leaf, cut	12328206	3
Caelo	Marshmallow root, nat...	270114	1
Caelo	Marshmallow root, pee...	12368906	1
Caelo	Marshmallow root, pee...	11101905	1
Caelo	Marshmallow root, pee...	15104503	1
Caelo	Marshmallow root, pee...	13057103	1
Caelo	Marshmallow root, pee...	12147311	7
Caelo	Marshmallow root, pee...	11101804	1
Caelo	Marshmallow root, pee...	26091201	1
Caelo	Marshmallow root, pee...	11101804	1
Caelo	Marshmallow root, pee...	12147311	1
Caelo	Marshmallow root, pee...	6071A130116	10
Caelo	Marshmallow root, pee...	11101907	7
Caelo	Masterwort root, cut	1210801	1
Caelo	Masterwort root, cut	13341201	1
Caelo	Masterwort root, cut	12101807	2
Caelo	Masterwort root, cut	13132002	1
Caelo	Masterwort root, cut	12101808	1
Caelo	Masterwort root, cut	12101806	2
Caelo	Maté leaf, cut	13061202	1
Caelo	Maté leaf, cut	12133801	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Maté leaf, cut	11091710	1
Caelo	Maté leaf, cut	12133803	3
Caelo	Maté leaf, cut	12133804	7
Caelo	Maté leaf, cut	12133807	7
Caelo	Maté leaf, cut	2141A120821	2
Caelo	Maté leaf, cut	13118003	4
Caelo	Maté leaf, cut	12133806	1
Caelo	Maté leaf, cut	13118001	1
Caelo	Matricaria flower	12244805	4
Caelo	Matricaria flower	12337608	1
Caelo	Matricaria flower	12244803	1
Caelo	Matricaria flower	19041306	1
Caelo	Matricaria flower	15258505	1
Caelo	Matricaria flower	12361701	3
Caelo	Matricaria flower	12244804	2
Caelo	Matricaria flower	12299701	1
Caelo	Matricaria flower	15258504	1
Caelo	Matricaria flower	22111202	1
Caelo	Matricaria flower	26081303	1
Caelo	Matricaria flower	13014501	1
Caelo	Matricaria flower	22111201	1
Caelo	Matricaria flower	12111903	2
Caelo	Matricaria flower	22011303	1
Caelo	Matricaria flower	12361702	6
Caelo	Matricaria flower	12299703	2
Caelo	Matricaria flower	13111201	1
Caelo	Matricaria flower	IntCh:313011	1
Caelo	Matricaria flower	15090003	2
Caelo	Matricaria flower	12294802	15
Caelo	Matricaria flower	13031306	1
Caelo	Matricaria flower	12337607	2
Caelo	Matricaria flower	16111202	1
Caelo	Matricaria flower	14276804	3
Caelo	Matricaria flower	19111207	1
Caelo	Matricaria flower	12337601	1
Caelo	Matricaria flower	12361707	2
Caelo	Matricaria flower	11101203	1
Caelo	Matricaria flower	12337610	1
Caelo	Matricaria flower	13210905	3
Caelo	Matricaria flower	11218501	1
Caelo	Matricaria flower	12361706	3
Caelo	Matricaria flower	12244803	3
Caelo	Matricaria flower	13210901-3478312171	1
Caelo	Matricaria flower	12294804	1
Caelo	Matricaria flower	12361704	1
Caelo	Matricaria flower	11169607	1
Caelo	Matricaria flower	13210903	1
Caelo	Matricaria flower	12361706	1
Caelo	Matricaria flower	13210901	3
Caelo	Matricaria flower	13210905	1
Caelo	Matricaria flower	13210906	1
Caelo	Matricaria flower	13011406	1
Caelo	Matricaria flower	13061201	1
Caelo	Matricaria flower	13288405	3
Caelo	Matricaria flower	17101204	1
Caelo	Matricaria flower	12361707	4
Caelo	Matricaria flower	10081201	1
Caelo	Matricaria flower	13210901-3546411013	1
Caelo	Matricaria flower	22081202	1
Caelo	Matricaria flower	12059606	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Matricaria flower	5081201	1
Caelo	Matricaria flower	20021313	1
Caelo	Matricaria flower	AKD5	1
Caelo	Matricaria flower	12244802	4
Caelo	Matricaria flower	27091203	1
Caelo	Matricaria flower	11214810	1
Caelo	Meadowsweet flower, r...	9326211	2
Caelo	Meadowsweet flower, r...	10349302	2
Caelo	Meadowsweet flower, r...	9326213	2
Caelo	Meadowsweet flower, r...	1227501	1
Caelo	Meadowsweet flower, r...	10349305	1
Caelo	Meadowsweet flower, r...	12275701	2
Caelo	Meadowsweet flower, r...	12275703	3
Caelo	Meadowsweet flower, r...	10349307	2
Caelo	Meadowsweet flower, r...	10349308	1
Caelo	Meadowsweet flower, r...	1225702	1
Caelo	Meadowsweet flower, r...	12275704-3372319032	1
Caelo	Meadowsweet flower, r...	12275704	4
Caelo	Meadowsweet flower, r...	10051301	1
Caelo	Meadowsweet flower, r...	12275702	1
Caelo	Meadowsweet herb, cut	12188703	1
Caelo	Meadowsweet herb, cut	12188702	1
Caelo	Meadowsweet herb, cut	62487517	1
Caelo	Meadowsweet herb, cut	12188803	10
Caelo	Meadowsweet herb, cut	12188801	1
Caelo	Meadowsweet herb, cut	4912E-01685	1
Caelo	Meadowsweet herb, cut	12188807	1
Caelo	Meadowsweet herb, cut	12188704	2
Caelo	Meadowsweet herb, cut	12188806	3
Caelo	Meadowsweet herb, cut	10381602	3
Caelo	Meadowsweet herb, cut	10381601	1
Caelo	Meadowsweet herb, cut	12188805	13
Caelo	Melatonin	163437/13	1
Caelo	Melatonin	18210311	2
Caelo	Melatonin	17305113	3
Caelo	Melatonin	17305118	2
Caelo	Melatonin	17305116	1
Caelo	Melatonin	17305114	3
Caelo	Melatonin	17/28-B04-345376	2
Caelo	Melatonin	17305119	3
Caelo	Melatonin	18210301	6
Caelo	Melatonin	17305122	4
Caelo	Melatonin	17305104	1
Caelo	Melatonin	18210304	5
Caelo	Melatonin	17350119	1
Caelo	Melatonin	18210306	4
Caelo	Melatonin	18210308	1
Caelo	Melatonin	17I28-B04-345376	1
Caelo	Melatonin	26051805	1
Caelo	Melatonin	17305105	7
Caelo	Melatonin	13167301	9
Caelo	Melatonin	16343711	2
Caelo	Melatonin	11267121	1
Caelo	Melatonin	14110305	1
Caelo	Melatonin	11267110	1
Caelo	Melatonin	14110308	3
Caelo	Melatonin	1126109	1
Caelo	Melatonin	11267122	1
Caelo	Melatonin	20021505	1
Caelo	Melatonin	13167314	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Melatonin	15093002	2
Caelo	Melatonin	15093001	6
Caelo	Melatonin	13167301	1
Caelo	Melatonin	17305110	5
Caelo	Melatonin	15093004	1
Caelo	Melatonin	14110314	2
Caelo	Melatonin	15093010	5
Caelo	Melatonin	12321415	4
Caelo	Melatonin	15093008	7
Caelo	Melatonin	15254306	7
Caelo	Melatonin	15093008	1
Caelo	Melatonin	15093013	2
Caelo	Melatonin	15254305	6
Caelo	Melatonin	15254305	1
Caelo	Melatonin	15093004	1
Caelo	Melatonin	15093012	2
Caelo	Melatonin	14110305	1
Caelo	Melatonin	14110302	1
Caelo	Melatonin	12321421	2
Caelo	Melatonin	120009	2
Caelo	Melatonin	12321411	1
Caelo	Melatonin	16121303	1
Caelo	Melatonin	12321416	5
Caelo	Melatonin	10071303	1
Caelo	Melatonin	11267126	2
Caelo	Melatonin	11267102	1
Caelo	Melatonin	12321422	4
Caelo	Melatonin	2313M03554	1
Caelo	Melatonin	11267109	1
Caelo	Melatonin	12321406	1
Caelo	Melatonin	20021312	1
Caelo	Melatonin	11267120	4
Caelo	Melatonin	13167310	1
Caelo	Melatonin	13167302	1
Caelo	Melatonin	13167309	1
Caelo	Melatonin	13167312	2
Caelo	Melatonin	11267119	2
Caelo	Melatonin	14110302	1
Caelo	Melatonin	1514E-03554	1
Caelo	Melatonin	13167312	1
Caelo	Melatonin	12321404	2
Caelo	Melatonin	15254307	15
Caelo	Melatonin	14110316	1
Caelo	Melatonin	150601	1
Caelo	Melatonin	15254327	1
Caelo	Melatonin	15254331	3
Caelo	Melatonin	15254328	2
Caelo	Melatonin	15254333	1
Caelo	Melatonin	15254332	2
Caelo	Melatonin	15254338	6
Caelo	Melatonin	15G06-B03-321298	5
Caelo	Melatonin	16266901	7
Caelo	Melatonin	16266908	4
Caelo	Melatonin	16266904	3
Caelo	Melatonin	15254308	8
Caelo	Melatonin	16266909	9
Caelo	Melatonin	16266903	1
Caelo	Melatonin	16266906	1
Caelo	Melatonin	16343706	1
Caelo	Melatonin	16343710	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Melatonin	16343707	3
Caelo	Melatonin	16343709	6
Caelo	Melatonin	16343705	1
Caelo	Melatonin	16343708	10
Caelo	Melatonin	17305101	2
Caelo	Melatonin	16343713	8
Caelo	Melatonin	17305102	14
Caelo	Melatonin	17305101-20180206Ma1	1
Caelo	Melatonin	17305106	6
Caelo	Melatonin	15254336	4
Caelo	Melatonin	15254324	7
Caelo	Melatonin	15254325	1
Caelo	Melatonin	15254313	2
Caelo	Melatonin	1524309	1
Caelo	Melatonin	160125M	1
Caelo	Melatonin	15254321	6
Caelo	Melatonin	15G06-B03-321299	3
Caelo	Melatonin	15254319	10
Caelo	Melatonin	15254320	3
Caelo	Melatonin	15254306	1
Caelo	Melilot herb, cut	25091301	1
Caelo	Melilot herb, cut	12289205	4
Caelo	Melilot herb, cut	11307606	6
Caelo	Melilot herb, cut	12289203	5
Caelo	Melilot herb, cut	19071201	1
Caelo	Melilot herb, cut	12289201	1
Caelo	Melilot herb, cut	22111303	1
Caelo	Melilot herb, cut	24081202	1
Caelo	Melilot herb, cut	12289206	1
Caelo	Melilot herb, cut	122891202	1
Caelo	Melilot herb, cut	8061204	1
Caelo	Melilot herb, cut	12231703	3
Caelo	Melilot herb, cut	12231701	3
Caelo	Melilot herb, cut	12289202	8
Caelo	Melilot herb, cut	107063	1
Caelo	Melilot herb, cut	12289207	1
Caelo	Melilot herb, cut	11307601	1
Caelo	Melilot herb, cut	11307607	1
Caelo	Melissa leaf, cut	12032303	1
Caelo	Melissa leaf, cut	13030803	2
Caelo	Melissa leaf, cut	10310919	1
Caelo	Melissa leaf, cut	19041303	1
Caelo	Melissa leaf, cut	12032307	2
Caelo	Melissa leaf, cut	12032310	1
Caelo	Melissa leaf, cut	13214720	1
Caelo	Melissa leaf, cut	13214707	1
Caelo	Melissa leaf, cut	13214705	1
Caelo	Melissa leaf, cut	24011407	1
Caelo	Melissa leaf, cut	12032319	5
Caelo	Melissa leaf, cut	14111319	1
Caelo	Melissa leaf, cut	12032316	3
Caelo	Melissa leaf, cut	12032309	1
Caelo	Melissa leaf, cut	25061303	1
Caelo	Melissa leaf, cut		1
Caelo	Melissa leaf, cut	CH.111111	1
Caelo	Melissa leaf, cut	13030805	1
Caelo	Melissa leaf, cut	1009	1
Caelo	Melissa leaf, cut	12032318	1
Caelo	Melissa leaf, cut	2161a121009	1
Caelo	Melissa leaf, cut	25021301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Menthol	14183012	3
Caelo	Menthol	15076115	2
Caelo	Menthol	15076104	1
Caelo	Menthol	515M-01855	1
Caelo	Menthol	1507601	1
Caelo	Menthol	15076215	1
Caelo	Menthol	15363805	9
Caelo	Menthol	15363804	9
Caelo	Menthol	15363803	2
Caelo	Menthol	15076112	10
Caelo	Menthol	15076106	10
Caelo	Menthol	14183004	13
Caelo	Menthol	814A-01855	1
Caelo	Menthol	15076113	8
Caelo	Menthol	14183020	1
Caelo	Menthol	14183006	1
Caelo	Menthol	14183014	2
Caelo	Menthol	15076101	8
Caelo	Menthol	14183017	3
Caelo	Menthol	12020411	10
Caelo	Menthol	14183011	1
Caelo	Menthol	14183008	3
Caelo	Menthol	13377412	4
Caelo	Menthol	13103802	1
Caelo	Menthol	13377413	1
Caelo	Menthol	Prüfmr. 3712U-01855	1
Caelo	Menthol	11202901	1
Caelo	Menthol	13377407	2
Caelo	Menthol	IntCh311052	1
Caelo	Menthol	750137	3
Caelo	Menthol	12246203	1
Caelo	Menthol	12246211	5
Caelo	Menthol	72173258	1
Caelo	Menthol	13159105	20
Caelo	Menthol	1320260	2
Caelo	Menthol	13159104	4
Caelo	Menthol	13103901	26
Caelo	Menthol	13G09-N09	2
Caelo	Menthol	314Q01855	1
Caelo	Menthol	29121103	1
Caelo	Menthol	10157615	1
Caelo	Menthol	25051301	1
Caelo	Menthol	7111201	1
Caelo	Menthol	687711	1
Caelo	Menthol	3912A-01855	1
Caelo	Menthol	1/AD1	1
Caelo	Menthol	25041402	1
Caelo	Menthol	13159110	2
Caelo	Menthol	19051401	1
Caelo	Menthol	13377404	15
Caelo	Menthol	13377402	3
Caelo	Menthol	13377408	2
Caelo	Menthol	12246205	8
Caelo	Menthol	1202041	1
Caelo	Menthol	11091405	1
Caelo	Menthol	13377409	9
Caelo	Menthol	10008029-02	1
Caelo	Menthol	20200712-2	1
Caelo	Menthol	12020410	1
Caelo	Menthol	130501001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Menthol	16095609	1
Caelo	Menthol	1046973	2
Caelo	Menthol	16081805	1
Caelo	Menthol	18152205	1
Caelo	Menthol	15286209	1
Caelo	Menthol	16095614	11
Caelo	Menthol	91(18128009)	1
Caelo	Menthol	18128009	5
Caelo	Menthol	18128007	1
Caelo	Menthol	1738804	1
Caelo	Menthol	170822M/16095616	1
Caelo	Menthol	171889035483	1
Caelo	Menthol	16095619	12
Caelo	Menthol	5117A-01855	1
Caelo	Menthol	17188906	2
Caelo	Menthol	16095604	1
Caelo	Menthol	17078006	1
Caelo	Menthol	160956217188904	1
Caelo	Menthol	90(17328804)	1
Caelo	Menthol	17328803	2
Caelo	Menthol	17351505	2
Caelo	Menthol	16095613	14
Caelo	Menthol	15286216	1
Caelo	Menthol	15076219	1
Caelo	Menthol	19091805	1
Caelo	Menthol	18128015	4
Caelo	Menthol	1608012-02	16
Caelo	Menthol	15286212	1
Caelo	Menthol	62470307	1
Caelo	Menthol	930716	1
Caelo	Menthol	160719M	1
Caelo	Menthol	15363810	12
Caelo	Menthol	16095602	3
Caelo	Menthol	16095620	6
Caelo	Menthol	15076118	1
Caelo	Menthol	17078004	1
Caelo	Menthol	15363815	12
Caelo	Menthol	1705022-01	27
Caelo	Menthol	92(18128015)	1
Caelo	Menthol	161202M/15076111	1
Caelo	Menthol	16095616	1
Caelo	Menthol	160956	1
Caelo	Menthol	16098614	1
Caelo	Menthol	2416	1
Caelo	Menthol	17188903	13
Caelo	Menthol	5315I-01855	1
Caelo	Menthol		1
Caelo	Menthol	1536815	1
Caelo	Menthol	16095603	4
Caelo	Menthol	16095606	13
Caelo	Mepivacaine hydrochlo. . .	13G05-N02	2
Caelo	Methacholine chloride	16329901	2
Caelo	Methionine	16274002	1
Caelo	Methionine	17257502	1
Caelo	Methoxsalen	17078108	3
Caelo	Methoxsalen	15244209	1
Caelo	Methoxsalen	17078101	2
Caelo	Methoxsalen	15244216	1
Caelo	Methoxsalen	15244211	1
Caelo	Methoxsalen	15244204	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Methoxsalen	11376433	1
Caelo	Methoxsalen	16/079/M	1
Caelo	Methoxsalen	15099802	1
Caelo	Methoxsalen	11376420	1
Caelo	Methoxsalen	1136428	1
Caelo	Methoxsalen	1409003-02	3
Caelo	Methoxsalen	1409003-01	9
Caelo	Methoxsalen	11376413	2
Caelo	Methyl-4-hydroxybenzo...	1251624	1
Caelo	Methyl-4-hydroxybenzo...	11345115	7
Caelo	Methyl-4-hydroxybenzo...	111003301	1
Caelo	Methyl-4-hydroxybenzo...	13268802	1
Caelo	Methyl-4-hydroxybenzo...	11090031	4
Caelo	Methyl-4-hydroxybenzo...	13268809	6
Caelo	Methyl-4-hydroxybenzo...	13268813	2
Caelo	Methyl-4-hydroxybenzo...	11090030	1
Caelo	Methyl-4-hydroxybenzo...	1108211930	1
Caelo	Methyl-4-hydroxybenzo...	9011303	1
Caelo	Methyl-4-hydroxybenzo...	13268804	6
Caelo	Methyl-4-hydroxybenzo...	713Q-01930	1
Caelo	Methyl-4-hydroxybenzo...	11345103	1
Caelo	Methyl-4-hydroxybenzo...	13081308	1
Caelo	Methyl-4-hydroxybenzo...	12251611	1
Caelo	Methyl-4-hydroxybenzo...	290114PML01	1
Caelo	Methyl-4-hydroxybenzo...	12251624	7
Caelo	Methyl-4-hydroxybenzo...	12251605	2
Caelo	Methyl-4-hydroxybenzo...	12251619	1
Caelo	Methyl-4-hydroxybenzo...	11090003	1
Caelo	Methyl-4-hydroxybenzo...	11345120	2
Caelo	Methyl-4-hydroxybenzo...	6091402	1
Caelo	Methyl-4-hydroxybenzo...	18176903	3
Caelo	Methyl-4-hydroxybenzo...	13268803	1
Caelo	Methyl-4-hydroxybenzo...	15063813	2
Caelo	Methyl-4-hydroxybenzo...	15298102	2
Caelo	Methyl-4-hydroxybenzo...	15063818	2
Caelo	Methyl-4-hydroxybenzo...	15063819	5
Caelo	Methyl-4-hydroxybenzo...	15063820	3
Caelo	Methyl-4-hydroxybenzo...	15444911	1
Caelo	Methyl-4-hydroxybenzo...	15063819/1210217	1
Caelo	Methyl-4-hydroxybenzo...	15444913	1
Caelo	Methyl-4-hydroxybenzo...	16256103	3
Caelo	Methyl-4-hydroxybenzo...	15444916	1
Caelo	Methyl-4-hydroxybenzo...	15444905	1
Caelo	Methyl-4-hydroxybenzo...	16256104	3
Caelo	Methyl-4-hydroxybenzo...	17158101	18
Caelo	Methyl-4-hydroxybenzo...	11090031	1
Caelo	Methyl-4-hydroxybenzo...	16256118	5
Caelo	Methyl-4-hydroxybenzo...	11345115	1
Caelo	Methyl-4-hydroxybenzo...	17343013	4
Caelo	Methyl-4-hydroxybenzo...	15298101	1
Caelo	Methyl-4-hydroxybenzo...	15298102	10
Caelo	Methyl-4-hydroxybenzo...	17343002	4
Caelo	Methyl-4-hydroxybenzo...	14120208	5
Caelo	Methyl-4-hydroxybenzo...	15063806	8
Caelo	Methyl-4-hydroxybenzo...	13268819	1
Caelo	Methyl-4-hydroxybenzo...	13268804	4
Caelo	Methyl-4-hydroxybenzo...	12081502	1
Caelo	Methyl-4-hydroxybenzo...	15063803	1
Caelo	Methyl-4-hydroxybenzo...	15063802	6
Caelo	Methyl-4-hydroxybenzo...	13268815	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Methyl-4-hydroxybenzo...	13268813	4
Caelo	Methyl-4-hydroxybenzo...	13268821	1
Caelo	Methyl-4-hydroxybenzo...	15222203	4
Caelo	Methylcellulose	15179602	11
Caelo	Methylcellulose	15179602	1
Caelo	Methylcellulose	14197502	2
Caelo	Methylcellulose	16060206	12
Caelo	Methylcellulose	14197514	5
Caelo	Methylcellulose	14197514	2
Caelo	Methylcellulose	16060210	10
Caelo	Methylcellulose	13130121	1
Caelo	Methylcellulose	15179605	7
Caelo	Methylcellulose	14197513	1
Caelo	Methylcellulose	16284006	6
Caelo	Methylcellulose	5081404	1
Caelo	Methylcellulose	16060202	1
Caelo	Methylcellulose	13130121	2
Caelo	Methylcellulose	13130115	5
Caelo	Methylcellulose	20021401	1
Caelo	Methylcellulose	13021401	1
Caelo	Methylcellulose	13130109	5
Caelo	Methylcellulose	12044931	2
Caelo	Methylcellulose	12044903	2
Caelo	Methylcellulose	10247604	1
Caelo	Methylcellulose	12044907	1
Caelo	Methylcellulose	12044918	2
Caelo	Methylcellulose	12044917	1
Caelo	Methylcellulose	14197502	4
Caelo	Methylcellulose	1628404	1
Caelo	Methylcellulose	13130104	5
Caelo	Methylcellulose	16284004	1
Caelo	Methylcellulose	17170405	5
Caelo	Methylcellulose	13130115	1
Caelo	Methylcellulose	18045101	3
Caelo	Methylcellulose	18045107	2
Caelo	Methylcellulose	18045114	2
Caelo	Methylcellulose	12044927	6
Caelo	Methylrosanilinium ch...	12H28-No6	1
Caelo	Methylrosanilinium ch...	9238106	1
Caelo	Methylrosanilinium ch...	9238107	1
Caelo	Methylrosanilinium ch...	13184801	1
Caelo	Methylrosanilinium ch...	12H28-N06	1
Caelo	Methylrosanilinium ch...	12/17-N05	1
Caelo	Metoclopramide hydroc...	15227001	15
Caelo	Metoclopramide hydroc...	3414U-018969	1
Caelo	Metoprolol tartrate	17G03-B09-351672	2
Caelo	Metronidazole	12339521	1
Caelo	Metronidazole	13095012	13
Caelo	Metronidazole	13095006	6
Caelo	Metronidazole	18071401	1
Caelo	Metronidazole	13095013	13
Caelo	Metronidazole	10368319	1
Caelo	Metronidazole	7071405	1
Caelo	Metronidazole	1210411	1
Caelo	Metronidazole	12339502	18
Caelo	Metronidazole	12110421	7
Caelo	Metronidazole	9715621	1
Caelo	Metronidazole	13095020	5
Caelo	Metronidazole	12110415	6
Caelo	Metronidazole	13214001	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Metronidazole	12110105	1
Caelo	Metronidazole	12091403	1
Caelo	Metronidazole	12339501	8
Caelo	Metronidazole	11319816	2
Caelo	Metronidazole	12110423	1
Caelo	Metronidazole	10368311	1
Caelo	Metronidazole	12339503	6
Caelo	Metronidazole	12339520	4
Caelo	Metronidazole	13095011	3
Caelo	Metronidazole	13095005	5
Caelo	Metronidazole	12L17-N11	19
Caelo	Metronidazole	11319825	2
Caelo	Metronidazole	10368310	1
Caelo	Metronidazole	12359501	1
Caelo	Metronidazole	12339519	13
Caelo	Metronidazole	12110414	3
Caelo	Metronidazole	12110405	8
Caelo	Metronidazole	12339511	15
Caelo	Metronidazole	28061204	1
Caelo	Metronidazole	12339509	4
Caelo	Metronidazole	3947/V04/16	1
Caelo	Metronidazole	13095001	10
Caelo	Metronidazole	11319824	2
Caelo	Metronidazole	607211870	1
Caelo	Metronidazole	1201156	1
Caelo	Metronidazole	11319819	1
Caelo	Metronidazole	11319807	1
Caelo	Metronidazole	2013U-01870	1
Caelo	Metronidazole	27092014A	1
Caelo	Metronidazole	121203001	3
Caelo	Metronidazole	12339508	4
Caelo	Metronidazole	20121121	4
Caelo	Metronidazole	121104085	1
Caelo	Metronidazole	2813E1870	1
Caelo	Metronidazole	11319810	2
Caelo	Metronidazole	10112907	1
Caelo	Metronidazole	8012014B	1
Caelo	Metronidazole	11h11-n03	1
Caelo	Metronidazole	12339507	6
Caelo	Metronidazole	1239504	1
Caelo	Metronidazole	102311870	1
Caelo	Metronidazole	12339504	6
Caelo	Metronidazole	12110403	1
Caelo	Metronidazole	20140612-12	1
Caelo	Metronidazole	2071304	1
Caelo	Metronidazole	15121202	1
Caelo	Metronidazole	1212030-21	1
Caelo	Metronidazole	7031401	1
Caelo	Metronidazole	13095004	11
Caelo	Metronidazole	12110411	5
Caelo	Metronidazole	13095002	10
Caelo	Metronidazole	614U1870	1
Caelo	Metronidazole	13095003	11
Caelo	Metronidazole	18051201	1
Caelo	Metronidazole	39/AKD16	1
Caelo	Metronidazole	10368313	1
Caelo	Metronidazole	12339526	1
Caelo	Metronidazole	12110413	4
Caelo	Metronidazole	3112M-01870	1
Caelo	Metronidazole	12110412	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Metronidazole	14111304	1
Caelo	Metronidazole	11319823	5
Caelo	Metronidazole	12339510	9
Caelo	Metronidazole	13214002	17
Caelo	Metronidazole	17290506	13
Caelo	Metronidazole	15215106	5
Caelo	Metronidazole	15215103	4
Caelo	Metronidazole	17290512	12
Caelo	Metronidazole	15A20-B07-303297	1
Caelo	Metronidazole	15215108	8
Caelo	Metronidazole	15215107	1
Caelo	Metronidazole	172905125723	1
Caelo	Metronidazole	15215104	5
Caelo	Metronidazole	12081507	1
Caelo	Metronidazole	1790512	1
Caelo	Metronidazole	13214029	8
Caelo	Metronidazole	13214021	1
Caelo	Metronidazole	13214025	7
Caelo	Metronidazole	13214018	5
Caelo	Metronidazole	14c14-b01-295545	1
Caelo	Metronidazole	180264001	1
Caelo	Metronidazole	17290507	7
Caelo	Metronidazole	17290509	3
Caelo	Metronidazole	15253001	11
Caelo	Metronidazole	15253003	8
Caelo	Metronidazole	15253015	10
Caelo	Metronidazole	1523006	1
Caelo	Metronidazole	15253011	9
Caelo	Metronidazole	15253006	7
Caelo	Metronidazole	15253005	3
Caelo	Metronidazole	1523013	1
Caelo	Metronidazole	15253002Ge22.12.15	1
Caelo	Metronidazole	17290504	13
Caelo	Metronidazole	15253004	13
Caelo	Metronidazole	1061804	2
Caelo	Metronidazole	15215105	1
Caelo	Metronidazole	18061803	1
Caelo	Metronidazole	17058013-20180619Ma6	1
Caelo	Metronidazole	17290503	12
Caelo	Metronidazole	15253013	9
Caelo	Metronidazole	17290515	3
Caelo	Metronidazole	1706022-01	8
Caelo	Metronidazole	14b27-b31-291230	1
Caelo	Metronidazole	17290513	2
Caelo	Metronidazole	13214024	3
Caelo	Metronidazole	113214021	1
Caelo	Metronidazole	13214015	8
Caelo	Metronidazole	13041503	1
Caelo	Metronidazole	13214017	7
Caelo	Metronidazole	13214016	2
Caelo	Metronidazole	13214012	4
Caelo	Metronidazole	13214011	1
Caelo	Metronidazole	13214013	1
Caelo	Metronidazole	13214014	9
Caelo	Metronidazole	132140022	1
Caelo	Metronidazole	2022015B	1
Caelo	Metronidazole	23011504	1
Caelo	Metronidazole	13095018	6
Caelo	Metronidazole	13214009	8
Caelo	Metronidazole	13095017	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Metronidazole	13214003	11
Caelo	Metronidazole	13214010	5
Caelo	Metronidazole	15253013/2016030911	1
Caelo	Metronidazole	1321409	1
Caelo	Metronidazole	19010203	2
Caelo	Metronidazole	14B27-B31-291230	1
Caelo	Metronidazole	39/AKD4	1
Caelo	Metronidazole	17290517	16
Caelo	Metronidazole	13214019	6
Caelo	Metronidazole	17290518	20
Caelo	Metronidazole	17290519	3
Caelo	Metronidazole	1790517	1
Caelo	Metronidazole	17290533	1
Caelo	Metronidazole	19011909	1
Caelo	Metronidazole	17290527	11
Caelo	Metronidazole	1729027	1
Caelo	Metronidazole	17290526	3
Caelo	Metronidazole	1739,527	1
Caelo	Metronidazole	49(17290526)	1
Caelo	Metronidazole	19010202	2
Caelo	Metronidazole	17290532	3
Caelo	Metronidazole	17290525	7
Caelo	Metronidazole	15253007	15
Caelo	Metronidazole	15253002	2
Caelo	Metronidazole	15253017	14
Caelo	Metronidazole	17058012	11
Caelo	Metronidazole	17058019	10
Caelo	Metronidazole	17058004	4
Caelo	Metronidazole	7091714	1
Caelo	Metronidazole	16089236	2
Caelo	Metronidazole	17058002	7
Caelo	Metronidazole	15253016	11
Caelo	Metronidazole	17058005,20171127Ma3	1
Caelo	Metronidazole	170580001	2
Caelo	Metronidazole	16089235	4
Caelo	Metronidazole	160892365302	1
Caelo	Metronidazole	16056235	1
Caelo	Metronidazole	16089228	11
Caelo	Metronidazole	16089229	10
Caelo	Metronidazole	16051703	1
Caelo	Metronidazole	16089222	7
Caelo	Metronidazole	17058001	8
Caelo	Metronidazole	16089123	1
Caelo	Metronidazole	17058013	15
Caelo	Metronidazole	17058012/2231117	1
Caelo	Metronidazole	280218010	1
Caelo	Metronidazole	1307002	1
Caelo	Metronidazole	17058013-20180404Ma4	1
Caelo	Metronidazole	170580165554	1
Caelo	Metronidazole	17058026	12
Caelo	Metronidazole	17058014	10
Caelo	Metronidazole	17058018	12
Caelo	Metronidazole	170580035444	1
Caelo	Metronidazole	17058011	7
Caelo	Metronidazole	170580123	1
Caelo	Metronidazole	1705801	1
Caelo	Metronidazole	1407003-2	1
Caelo	Metronidazole	17058015	10
Caelo	Metronidazole	17058027	7
Caelo	Metronidazole	17058010	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Metronidazole	16089231	2
Caelo	Metronidazole	17058017	9
Caelo	Metronidazole	16089223	4
Caelo	Metronidazole	16089228-20170818Ma4	1
Caelo	Metronidazole	16089224	8
Caelo	Metronidazole	15253028	9
Caelo	Metronidazole	15253027	11
Caelo	Metronidazole	15253029	3
Caelo	Metronidazole	15253026	19
Caelo	Metronidazole	1523026	1
Caelo	Metronidazole	15253024	11
Caelo	Metronidazole	15253022	4
Caelo	Metronidazole	ch.b.15253027	1
Caelo	Metronidazole	15258016	1
Caelo	Metronidazole	16089220	14
Caelo	Metronidazole	16253015	1
Caelo	Metronidazole	15253020	5
Caelo	Metronidazole	13095019	7
Caelo	Metronidazole	1523017	1
Caelo	Metronidazole	1525016	1
Caelo	Metronidazole	15253018	6
Caelo	Metronidazole	15253021	5
Caelo	Metronidazole	16089202	9
Caelo	Metronidazole	15253019	3
Caelo	Metronidazole	15253033	5
Caelo	Metronidazole	16089209	13
Caelo	Metronidazole	16089221	4
Caelo	Metronidazole	1619-b04-326150	1
Caelo	Metronidazole	16089225	1
Caelo	Metronidazole	16089219	5
Caelo	Metronidazole	23021713	1
Caelo	Metronidazole	16089215	12
Caelo	Metronidazole	16089217	8
Caelo	Metronidazole	16089215/1140217	1
Caelo	Metronidazole	16089216	6
Caelo	Metronidazole	10112903	1
Caelo	Metronidazole	16089213	14
Caelo	Metronidazole	15253032,20161215Ma8	1
Caelo	Metronidazole	160902801	2
Caelo	Metronidazole	16/19-B04-326155	3
Caelo	Metronidazole	16089218	8
Caelo	Metronidazole	16089212	16
Caelo	Miconazole nitrate	12k02n06	1
Caelo	Miconazole nitrate	11324913	1
Caelo	Miconazole nitrate	12243805	7
Caelo	Miconazole nitrate	11104108	1
Caelo	Miconazole nitrate	12243806	1
Caelo	Miconazole nitrate	1/AD1	1
Caelo	Miconazole nitrate		1
Caelo	Miconazole nitrate	12374308	6
Caelo	Miconazole nitrate	9177904	1
Caelo	Miconazole nitrate	28051306	1
Caelo	Miconazole nitrate	12374307	16
Caelo	Miconazole nitrate	12243802	5
Caelo	Miconazole nitrate	11324901	1
Caelo	Miconazole nitrate	11324909	1
Caelo	Miconazole nitrate	12243801	4
Caelo	Miconazole nitrate	12243803	11
Caelo	Miconazole nitrate	11324912	2
Caelo	Miconazole nitrate	11324905	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Miconazole nitrate	11104103	2
Caelo	Miconazole nitrate	11104114	1
Caelo	Miconazole nitrate	11324902	1
Caelo	Miconazole nitrate	11104109	2
Caelo	Miconazole nitrate	11104113	1
Caelo	Miconazole nitrate	4061304	1
Caelo	Miconazole nitrate	11243803	1
Caelo	Miconazole nitrate	13185603	11
Caelo	Miconazole nitrate	13352811	6
Caelo	Miconazole nitrate	13185602	7
Caelo	Miconazole nitrate	14149311	16
Caelo	Miconazole nitrate	13352811Mg:5g	1
Caelo	Miconazole nitrate	16144704	14
Caelo	Miconazole nitrate	1419310	1
Caelo	Miconazole nitrate	16144710	13
Caelo	Miconazole nitrate	16144717	8
Caelo	Miconazole nitrate	16144716	25
Caelo	Miconazole nitrate	16144711	9
Caelo	Miconazole nitrate	16e04b03327367	1
Caelo	Miconazole nitrate	16144721	11
Caelo	Miconazole nitrate	16E04-B03-325448	1
Caelo	Miconazole nitrate	29001614	1
Caelo	Miconazole nitrate	17217505	16
Caelo	Miconazole nitrate	17217507	1
Caelo	Miconazole nitrate	17C01B02340778	1
Caelo	Miconazole nitrate	18057904	1
Caelo	Miconazole nitrate	14149308	13
Caelo	Miconazole nitrate	122438	1
Caelo	Miconazole nitrate	17217503	6
Caelo	Miconazole nitrate	18057903	10
Caelo	Miconazole nitrate	18057909	1
Caelo	Miconazole nitrate	18057906	12
Caelo	Miconazole nitrate	18057902	8
Caelo	Miconazole nitrate	16144702	22
Caelo	Miconazole nitrate	16e04-b03-321367	1
Caelo	Miconazole nitrate	14149309	7
Caelo	Miconazole nitrate	13185601	4
Caelo	Miconazole nitrate	15137902	17
Caelo	Miconazole nitrate	15137906	1
Caelo	Miconazole nitrate	15137903	18
Caelo	Miconazole nitrate	15137911	1
Caelo	Miconazole nitrate	15I08-B07-312317	3
Caelo	Miconazole nitrate	151379	1
Caelo	Miconazole nitrate	15137907	9
Caelo	Miconazole nitrate	15137905	12
Caelo	Miconazole nitrate	2900/11/15	3
Caelo	Miconazole nitrate	15137908	7
Caelo	Miconazole nitrate	15A13-B07	8
Caelo	Miconazole nitrate	18057916	1
Caelo	Miconazole nitrate	15I08-B07-315481	2
Caelo	Miconazole nitrate	15137909	9
Caelo	Miconazole nitrate	14149318	19
Caelo	Miconazole nitrate	15/05-B07-312317	1
Caelo	Miconazole nitrate	15137910	7
Caelo	Miconazole nitrate	15137913	4
Caelo	Miconazole nitrate	15137914	4
Caelo	Miconazole nitrate	15137915	5
Caelo	Miconazole nitrate	15108-B07-314838	1
Caelo	Miconazole nitrate	15108-B07-318686	1
Caelo	Miconazole nitrate	16144701	27

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Miconazole nitrate	16144703	17
Caelo	Miconazole nitrate	2900/10/14	2
Caelo	Miconazole nitrate	18057912	9
Caelo	Miconazole nitrate	16144726	15
Caelo	Miconazole nitrate	18B06_b03_34778	1
Caelo	Miconazole nitrate	14149310	8
Caelo	Miconazole nitrate	14149303	5
Caelo	Miconazole nitrate	133528Mg:5g	1
Caelo	Miconazole nitrate	1122014A	1
Caelo	Miconazole nitrate	14149301	17
Caelo	Miconazole nitrate	1419301	1
Caelo	Miconazole nitrate	14149302	15
Caelo	Miconazole nitrate	14149305	1
Caelo	Miconazole nitrate	14B18-B02-291233	1
Caelo	Miconazole nitrate	14D25-B06-294960	3
Caelo	Miconazole nitrate	13352806	12
Caelo	Miconazole nitrate	13352802	21
Caelo	Miconazole nitrate	13342801	1
Caelo	Miconazole nitrate	13352805	1
Caelo	Miconazole nitrate	411092-2-P	1
Caelo	Miconazole nitrate	13352801	21
Caelo	Miconazole nitrate	13185610	6
Caelo	Miconazole nitrate	13185609	9
Caelo	Miconazole nitrate	13185611	18
Caelo	Miconazole nitrate	AKD9	1
Caelo	Miconazole nitrate	13185612	5
Caelo	Miconazole nitrate	18057907	6
Caelo	Miconazole nitrate	14149317	9
Caelo	Miconazole nitrate	18057914	3
Caelo	Miconazole nitrate	18057913	4
Caelo	Miconazole nitrate	18282105	1
Caelo	Miconazole nitrate	18282106	1
Caelo	Miconazole nitrate	18282109	1
Caelo	Minoxidil	13179504	2
Caelo	Minoxidil	13373106	1
Caelo	Minoxidil	14234505	6
Caelo	Minoxidil	14234502	2
Caelo	Minoxidil	15436414	4
Caelo	Minoxidil	14062018-RQ-A	1
Caelo	Minoxidil	170500801	9
Caelo	Minoxidil	14093706	3
Caelo	Minoxidil	15230804	4
Caelo	Minoxidil	15230812	1
Caelo	Minoxidil	16J24-B05-331128	1
Caelo	Minoxidil	15230802	2
Caelo	Minoxidil	22031402	1
Caelo	Minoxidil	15230802	1
Caelo	Minoxidil	1807018-01	21
Caelo	Minoxidil	17011602	1
Caelo	Minoxidil	1203035	2
Caelo	Minoxidil	12H07-N03	15
Caelo	Minoxidil	13179501-3467312121	1
Caelo	Minoxidil	13/23-N07	1
Caelo	Minoxidil	20160211-1	1
Caelo	Minoxidil	15230805	2
Caelo	Minoxidil	15436412	1
Caelo	Minoxidil	15E19-B06-309180	1
Caelo	Minoxidil	17H08-B06-342282	1
Caelo	Minoxidil	12273710	1
Caelo	Minoxidil	14L03-B02-301431	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Minoxidil	13179503	6
Caelo	Minoxidil	15E19-B06-308183	1
Caelo	Minoxidil	14093707	1
Caelo	Minoxidil	13179501	4
Caelo	Minoxidil	150401301	1
Caelo	Minoxidil	3013M-02844	1
Caelo	Minoxidil	18016107	2
Caelo	Minoxidil	13039503	3
Caelo	Minoxidil	1309501	1
Caelo	Minoxidil	1070318/17227803	1
Caelo	Minoxidil	K0071/16/04/065	1
Caelo	Minoxidil	180701801	1
Caelo	Minoxidil	1901005-01	10
Caelo	Minoxidil	18149802	3
Caelo	Minoxidil	2905047-01	1
Caelo	Minoxidil	14093703	4
Caelo	Minoxidil	13039501	8
Caelo	Minoxidil	12273703	2
Caelo	Minoxidil	17227801	9
Caelo	Minoxidil	17227803	4
Caelo	Minoxidil	12273704	3
Caelo	Minoxidil	5081502	1
Caelo	Minoxidil	16J24-B05-332092	4
Caelo	Minoxidil	5217M-02844	1
Caelo	Minoxidil	14234509	3
Caelo	Minoxidil	13K04B05290380	1
Caelo	Minoxidil	12023908	1
Caelo	Minoxidil	170043-01	1
Caelo	Minoxidil	apk007.1-16-04-065-1	1
Caelo	Minoxidil	31051305	1
Caelo	Minoxidil	15436415	1
Caelo	Minoxidil	1411006-01	41
Caelo	Minoxidil	13179507	3
Caelo	Minoxidil	121005601	1
Caelo	Minoxidil	ch.b.15217502	1
Caelo	Minoxidil	ch.b.15436204	1
Caelo	Minoxidil	ch.b.15436202	1
Caelo	Minoxidil	1509035	1
Caelo	Minoxidil	12023901	2
Caelo	Minoxidil	15230811	3
Caelo	Minoxidil	15230811	1
Caelo	Minoxidil	11302003	1
Caelo	Minoxidil	15436204	2
Caelo	Minoxidil	160013-01	1
Caelo	Minoxidil	15436202	5
Caelo	Minoxidil	160707-01	1
Caelo	Minoxidil	160201302	2
Caelo	Minoxidil	15K17-B05-318692	1
Caelo	Minoxidil	15230803	3
Caelo	Minoxidil	AP011/16M	2
Caelo	Minoxidil	11302005	1
Caelo	Minoxidil	15436202/1060217	1
Caelo	Minoxidil	1520002455	1
Caelo	Minoxidil	12023902	1
Caelo	Minoxidil	11150207	1
Caelo	Minoxidil	15436201	5
Caelo	Minoxidil	15436404	4
Caelo	Minoxidil	20151112-5	1
Caelo	Minoxidil	15K17-B06-321926	7
Caelo	Minoxidil	15k17-b06-324965	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Minoxidil	11K11-n09	1
Caelo	Minoxidil	15436202	1
Caelo	Minoxidil	11e17-n07	1
Caelo	Minoxidil	15230806	1
Caelo	Minoxidil	12023907	4
Caelo	Minoxidil	15H07-B03-317340	1
Caelo	Minoxidil	11150209	1
Caelo	Mistletoe, cut	29071302	1
Caelo	Mistletoe, cut	11294004	5
Caelo	Mistletoe, cut	11266712	8
Caelo	Mistletoe, cut	29081302	1
Caelo	Mistletoe, cut	10273420	5
Caelo	Mistletoe, cut	11266714	4
Caelo	Mistletoe, cut	3813E-01695	1
Caelo	Mistletoe, cut	11266709	9
Caelo	Mistletoe, cut	29061201	1
Caelo	Mistletoe, cut	10273418	1
Caelo	Mistletoe, cut	17061304	1
Caelo	Mistletoe, cut	20121303	1
Caelo	Mistletoe, cut	14292506	2
Caelo	Mistletoe, cut	10081205	1
Caelo	Mistletoe, cut	27121305	1
Caelo	Mistletoe, cut	6061306	1
Caelo	Mistletoe, cut	10273412	1
Caelo	Mistletoe, cut	11266702	4
Caelo	Mistletoe, cut	11266704	14
Caelo	Mistletoe, cut	11266703	3
Caelo	Mistletoe, cut	12011301	1
Caelo	Mistletoe, cut	26071302	1
Caelo	Mistletoe, cut	11266710	7
Caelo	Mistletoe, cut	11294005	15
Caelo	Mistletoe, cut	11266708	8
Caelo	Mistletoe, cut	23111204	1
Caelo	Mistletoe, cut	24011401	1
Caelo	Mistletoe, cut	11021303	1
Caelo	Mistletoe, cut	11294003	2
Caelo	Mistletoe, cut	11266713	2
Caelo	Mistletoe, cut	11266712-3450311162	1
Caelo	Mistletoe, cut	29101302	1
Caelo	Mistletoe, cut	4313E-01695	1
Caelo	Mistletoe, cut	10091201	1
Caelo	Mistletoe, cut	14292514	1
Caelo	Mistletoe, cut	11266705	3
Caelo	Mistletoe, cut	28101304	1
Caelo	Mistletoe, cut	18031303	1
Caelo	Mistletoe, cut	10273407	1
Caelo	Molecular distilled g...	11160012	1
Caelo	Mometasone furoate, m...	14C07-B05-291268	12
Caelo	Mometasone furoate, m...	16028614	1
Caelo	Mometasone furoate, m...	18026006	2
Caelo	Mometasone furoate, m...	18026003	2
Caelo	Mometasone furoate, m...	1635206	1
Caelo	Mometasone furoate, m...	13J23-N02	17
Caelo	Mometasone furoate, m...	16353204	2
Caelo	Mometasone furoate, m...	NT1568M	4
Caelo	Mometasone furoate, m...	18B19B07-351001	1
Caelo	Mometasone furoate, m...	14E20-B04	3
Caelo	Mometasone furoate, m...	12351207	1
Caelo	Mometasone furoate, m...	16353206	10
Caelo	Mometasone furoate, m...	18026009	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Mometasone furoate, m...	16353210	1
Caelo	Mometasone furoate, m...	16353205	2
Caelo	Mometasone furoate, m...	18334104	2
Caelo	Mometasone furoate, m...	18026011	8
Caelo	Mometasone furoate, m...	1635209	1
Caelo	Mometasone furoate, m...	701834-0002	1
Caelo	Mometasone furoate, m...	17348901	3
Caelo	Mometasone furoate, m...	18026007	1
Caelo	Mometasone furoate, m...	18026001	12
Caelo	Mometasone furoate, m...	16353209	13
Caelo	Mometasone furoate, m...	13304604	1
Caelo	Mometasone furoate, m...	12339604	1
Caelo	Mometasone furoate, m...	17D19-B02-34476	1
Caelo	Mometasone furoate, m...	16353208	7
Caelo	Mometasone furoate, m...	12351202	1
Caelo	Mometasone furoate, m...	12351215	2
Caelo	Mometasone furoate, m...	16353207	11
Caelo	Mometasone furoate, m...	18026004	6
Caelo	Mometasone furoate, m...	16353211	8
Caelo	Mometasone furoate, m...	18026017	2
Caelo	Mometasone furoate, m...	12035601	2
Caelo	Mometasone furoate, m...	14281102	3
Caelo	Mometasone furoate, m...	NV1000M	1
Caelo	Mometasone furoate, m...	15189804	1
Caelo	Mometasone furoate, m...	15/11-Bo2-314479	1
Caelo	Mometasone furoate, m...	15189809	2
Caelo	Mometasone furoate, m...	16181715	7
Caelo	Mometasone furoate, m...	16028606	5
Caelo	Mometasone furoate, m...	12035607	2
Caelo	Mometasone furoate, m...	240913a	1
Caelo	Mometasone furoate, m...	16028601	2
Caelo	Mometasone furoate, m...	1510036-001	1
Caelo	Mometasone furoate, m...	13F27-N02	7
Caelo	Mometasone furoate, m...	12035608	1
Caelo	Mometasone furoate, m...	15189807	3
Caelo	Mometasone furoate, m...	12351205	1
Caelo	Mometasone furoate, m...	15189802	1
Caelo	Mometasone furoate, m...	15110305	3
Caelo	Mometasone furoate, m...	15b11b03	1
Caelo	Mometasone furoate, m...	14281103	1
Caelo	Mometasone furoate, m...	14281101	2
Caelo	Mometasone furoate, m...	15279004	4
Caelo	Mometasone furoate, m...	16028609	6
Caelo	Mometasone furoate, m...	200315A	1
Caelo	Mometasone furoate, m...	16181705	4
Caelo	Mometasone furoate, m...	16181710	1
Caelo	Mometasone furoate, m...	16181704	4
Caelo	Mometasone furoate, m...	13304605	1
Caelo	Mometasone furoate, m...	18026015	2
Caelo	Mometasone furoate, m...	16181709	1
Caelo	Mometasone furoate, m...	16028605	4
Caelo	Mometasone furoate, m...	16A20B06319748	1
Caelo	Mometasone furoate, m...	16181708	3
Caelo	Mometasone furoate, m...	18102013B	1
Caelo	Mometasone furoate, m...	16C24_B03_324852	2
Caelo	Mometasone furoate, m...	16028602	5
Caelo	Mometasone furoate, m...	15G08-B11-322739	1
Caelo	Mometasone furoate, m...	1704012A	1
Caelo	Mometasone furoate, m...	16028603	2
Caelo	Mometasone furoate, m...	14H06-B07298481	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Mometasone furoate, m...	N01352M	1
Caelo	Mometasone furoate, m...	16028607	6
Caelo	Mometasone furoate, m...	16181714	1
Caelo	Mometasone furoate, m...	13304606	1
Caelo	Mometasone furoate, m...	16J24-B11-330978	5
Caelo	Mometasone furoate, m...		1
Caelo	Motherwort herb, cut	12376408	2
Caelo	Motherwort herb, cut	12376403	1
Caelo	Motherwort herb, cut	12376407	1
Caelo	Motherwort herb, cut	11293004	1
Caelo	Motherwort herb, cut	5061203	1
Caelo	Mullein flower, cut	11020714	4
Caelo	Mullein flower, cut	11233401	4
Caelo	Mullein flower, cut	11366402	1
Caelo	Mullein flower, cut	11233402	1
Caelo	Mullein flower, cut	12192901	3
Caelo	Mullein flower, cut	12192903	3
Caelo	Mullein flower, cut	18031301	1
Caelo	Mullein flower, cut	16081201	1
Caelo	Mullein flower, cut	11333402	1
Caelo	Mullein flower, cut	11366403	1
Caelo	Mullein flower, cut	12192904	1
Caelo	Mullein flower, cut	12192907	3
Caelo	Mullein flower, cut	12170901	1
Caelo	Mullein flower, cut	1219203	1
Caelo	Mullein flower, cut	12317104	1
Caelo	Mullein flower, cut	12170906	1
Caelo	Mullein flower, cut	521505390	1
Caelo	Mullein flower, cut	13229702	1
Caelo	Mullein flower, cut	12192905	5
Caelo	Mullein flower, cut	11020713	2
Caelo	Mullein flower, cut	13229701	2
Caelo	Mullein flower, cut	11020902	1
Caelo	Mullein flower, cut	11020802	2
Caelo	Mupirocin calcium	17356502	1
Caelo	Mupirocin calcium	24041801	1
Caelo	Mupirocin calcium	83200301013	1
Caelo	Mupirocin calcium	9041801	1
Caelo	Myrrh	11276806	1
Caelo	Myrrh	11277111	1
Caelo	Myrrh	11277105	1
Caelo	Myrrh	11277111	1
Caelo	Myrrh	12244404	1
Caelo	Myrrh	12269904	4
Caelo	Myrrh	14347509	1
Caelo	Myrrh	30121413	1
Caelo	Myrrh	28101401	1
Caelo	Myrrh	25071404	1
Caelo	Myrrh	14244602	1
Caelo	Myrrh	12269907	1
Caelo	Myrrh	12269913	1
Caelo	Myrrh	13108701	2
Caelo	Myrrh	15312606	1
Caelo	Myrrh	12244406	1
Caelo	Myrrh	12269901	2
Caelo	Myrrh	11276810	1
Caelo	Myrrh	15021410	1
Caelo	Naphazoline hydrochlo...	1307020-02	1
Caelo	Naphazoline hydrochlo...	15223301	2
Caelo	Naphazoline hydrochlo...	15071404	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Naphazoline hydrochlo...	12375203	1
Caelo	Naphazoline hydrochlo...	1216203	1
Caelo	Naphazoline hydrochlo...	12176203	2
Caelo	Naphazoline hydrochlo...	18119302	1
Caelo	Naphazoline hydrochlo...	2275310	1
Caelo	Naphazoline hydrochlo...	732570001	1
Caelo	Naphazoline hydrochlo...	12375211	1
Caelo	Naphazoline hydrochlo...	15J05-B05/171025NHC	1
Caelo	Naphazoline hydrochlo...	15223303	1
Caelo	Naphazoline hydrochlo...	11081107	1
Caelo	Naphazoline hydrochlo...	16174101	3
Caelo	Naphazoline hydrochlo...	1506015-2	1
Caelo	Naphazoline hydrochlo...	11081107	1
Caelo	Naphazoline hydrochlo...	9285103	3
Caelo	Naphazoline hydrochlo...	140703402	1
Caelo	Naphazoline hydrochlo...	1206037	1
Caelo	Neomycin sulfate	5/AKD35	1
Caelo	Neomycin sulfate	563313062905	1
Caelo	Neomycin sulfate	13D10N05	1
Caelo	Neomycin sulfate	5/AKD36	1
Caelo	Neomycin sulfate	24111203	1
Caelo	Neomycin sulfate	15298902	3
Caelo	Neomycin sulfate	13062908	1
Caelo	Neomycin sulfate	12276410	2
Caelo	Neomycin sulfate	13062923	3
Caelo	Neomycin sulfate	13G09-N07	14
Caelo	Neomycin sulfate	13062922	1
Caelo	Neomycin sulfate	13062923	2
Caelo	Neomycin sulfate	15237602	6
Caelo	Neomycin sulfate	12276405	1
Caelo	Neomycin sulfate	130303303	3
Caelo	Neomycin sulfate	9051402	1
Caelo	Neomycin sulfate	18161303	3
Caelo	Neomycin sulfate	16165008	1
Caelo	Neomycin sulfate	16165016	9
Caelo	Neomycin sulfate	9105316	1
Caelo	Neomycin sulfate	18A25-B05-346135	1
Caelo	Neomycin sulfate	16165004	1
Caelo	Neomycin sulfate	12175401	4
Caelo	Neomycin sulfate	16B10-B50-328165	7
Caelo	Neomycin sulfate	16122612	5
Caelo	Neomycin sulfate	16122611	2
Caelo	Neomycin sulfate	10157505	5
Caelo	Neomycin sulfate	13062921	1
Caelo	Neomycin sulfate	1601046	1
Caelo	Neomycin sulfate	16122601	1
Caelo	Neomycin sulfate	16165018	4
Caelo	Neomycin sulfate	16165015	1
Caelo	Neomycin sulfate	16122618	1
Caelo	Neomycin sulfate	161650145351	1
Caelo	Neostigmine bromide	15228103	1
Caelo	Neostigmine bromide	15228101	1
Caelo	Neostigmine bromide	13177201	2
Caelo	Neostigmine bromide	13177206	5
Caelo	Neostigmine bromide	13177203	2
Caelo	Neostigmine bromide	131772303B	1
Caelo	Neostigmine bromide	12172304	1
Caelo	Neostigmine bromide	14256802	3
Caelo	Neostigmine bromide	131772303A	1
Caelo	Nettle leaf, cut	12034403	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Nettle leaf, cut	10375901	2
Caelo	Nettle leaf, cut	12183908	5
Caelo	Nettle leaf, cut	20111306	1
Caelo	Nettle leaf, cut	12184702	2
Caelo	Nettle leaf, cut	12183908-3443311112	1
Caelo	Nettle leaf, cut	33/AD1	1
Caelo	Nettle leaf, cut	12184706	3
Caelo	Nettle leaf, cut	intch:313011	1
Caelo	Nettle leaf, cut	12183910	11
Caelo	Nettle leaf, cut	1283910	1
Caelo	Nettle leaf, cut	29071303	1
Caelo	Nettle leaf, cut	12184005	7
Caelo	Nettle leaf, cut	72085497	1
Caelo	Nettle leaf, cut	83374109	1
Caelo	Nettle leaf, cut	12814001	1
Caelo	Nettle leaf, cut	7031302	1
Caelo	Nettle leaf, cut	83142389	1
Caelo	Nettle leaf, cut	2051302	1
Caelo	Nettle leaf, cut	12183908-3391311070	1
Caelo	Nettle leaf, cut	2381A130122	24
Caelo	Nettle leaf, cut	12184004	8
Caelo	Nettle leaf, cut	12183911	16
Caelo	Nettle leaf, cut	12183909	1
Caelo	Nettle leaf, cut	12183908-3405310181	1
Caelo	Nettle leaf, cut	12184001	6
Caelo	Nettle leaf, cut	13091302	1
Caelo	Nettle leaf, cut	12183914	8
Caelo	Nettle leaf, cut	10291422	3
Caelo	Nettle leaf, cut	12184502	2
Caelo	Nettle leaf, cut	12183901	2
Caelo	Nettle leaf, cut	10291405	2
Caelo	Nettle leaf, cut	10282807	2
Caelo	Nettle leaf, cut	14344907	1
Caelo	Nettle leaf, cut	12195002	2
Caelo	Nettle leaf, cut	13386801	1
Caelo	Nettle leaf, cut	10282806	1
Caelo	Nettle leaf, cut	10291414	3
Caelo	Nettle leaf, cut	12034414	1
Caelo	Nettle leaf, cut		3
Caelo	Nettle leaf, cut	18071202	1
Caelo	Nettle leaf, cut	10291413	1
Caelo	Nettle leaf, cut	120375901	1
Caelo	Nettle leaf, cut	10162712	1
Caelo	Nettle leaf, cut	12034402	3
Caelo	Nettle leaf, cut	3021401	1
Caelo	Nettle leaf, cut	17121205	1
Caelo	Nettle leaf, cut	12195101	3
Caelo	Nettle leaf, cut	10291411	2
Caelo	Nettle leaf, cut	12184503	3
Caelo	Nettle root, cut	6441A121024	1
Caelo	Nettle root, cut	11254904	3
Caelo	Nettle root, cut	11265004	4
Caelo	Nettle root, cut	6631B130610	1
Caelo	Nettle root, cut	11254903	5
Caelo	Nettle root, cut	11254902	3
Caelo	Nettle root, cut	11265006	1
Caelo	Nettle root, cut	10081202	1
Caelo	Nettle root, cut	11265002	5
Caelo	Nettle root, cut	11265003	2
Caelo	Nettle seed	16327003	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Nettle seed	20111301	1
Caelo	Nettle seed	13010612	2
Caelo	Nettle seed	15290602	1
Caelo	Nettle seed	16327001	2
Caelo	Nettle seed	12131505	1
Caelo	Nettle seed	13010607	1
Caelo	Nettle seed	17282601	1
Caelo	Nicotinamide	1511010-01	5
Caelo	Nicotinamide	13097504	1
Caelo	Nicotinamide	13097506	1
Caelo	Nicotinamide	12268302	1
Caelo	Nicotinamide	13352002	1
Caelo	Nicotinamide	14246806	1
Caelo	Nicotinamide	17104108	1
Caelo	Nicotinamide	16304803	1
Caelo	Non-ionic emulsifying...	12254402	1
Caelo	Non-ionic emulsifying...	311625	3
Caelo	Non-ionic emulsifying...	21091302	1
Caelo	Non-ionic emulsifying...	15287902	3
Caelo	Non-ionic emulsifying...	14241503	1
Caelo	Norethisterone acetate	15417001	1
Caelo	Norethisterone acetate	15417008	1
Caelo	Nutmeg	12376802	2
Caelo	Nutmeg	12376805	1
Caelo	Nutmeg, powdered (spi...	12102305	2
Caelo	Nystatin	12017010	1
Caelo	Nystatin	13101202	1
Caelo	Nystatin	12158021	2
Caelo	Nystatin	18072112	2
Caelo	Nystatin	18072014	1
Caelo	Nystatin	18072114	6
Caelo	Nystatin	18072108	2
Caelo	Nystatin	18072109	3
Caelo	Nystatin	1612001-03	1
Caelo	Nystatin	18072105	3
Caelo	Nystatin	11069706	2
Caelo	Nystatin	18072113	7
Caelo	Nystatin	12158005	1
Caelo	Nystatin	6082013B	1
Caelo	Nystatin	14363514	1
Caelo	Nystatin	13101201	1
Caelo	Nystatin	12158009	1
Caelo	Nystatin	2091304	1
Caelo	Nystatin	121580921	1
Caelo	Nystatin	12158022	1
Caelo	Nystatin	12138021	2
Caelo	Nystatin	12158020	1
Caelo	Nystatin	15335310	1
Caelo	Nystatin	12158019	3
Caelo	Nystatin	12158010	4
Caelo	Nystatin	11069731	1
Caelo	Nystatin	11069704	1
Caelo	Nystatin	12017003	1
Caelo	Nystatin	3101201	1
Caelo	Nystatin	12017011	2
Caelo	Nystatin	12158006	7
Caelo	Nystatin	12158004	3
Caelo	Nystatin	C129330	1
Caelo	Nystatin	12017009	1
Caelo	Nystatin	11069730	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Nystatin	12017004	3
Caelo	Nystatin	12158004	1
Caelo	Nystatin	19028801	2
Caelo	Nystatin	1 (18072128)	1
Caelo	Nystatin	18072124	4
Caelo	Nystatin	18072120	2
Caelo	Nystatin	11069716	1
Caelo	Nystatin	18072123	6
Caelo	Nystatin	12017015	2
Caelo	Nystatin	11069719	1
Caelo	Nystatin	13173401	11
Caelo	Oak bark, cut	11351209	9
Caelo	Oak bark, cut	10348613	2
Caelo	Oak bark, cut	11351208	3
Caelo	Oak bark, cut	14111316	1
Caelo	Oak bark, cut	10348605	1
Caelo	Oak bark, cut	11351203	1
Caelo	Oak bark, cut	11351207-3384310110	1
Caelo	Oak bark, cut	123test	1
Caelo	Oak bark, cut	161212	1
Caelo	Oak bark, cut	15146110	1
Caelo	Oak bark, cut	15146105	1
Caelo	Oak bark, cut	15146105	4
Caelo	Oak bark, cut	15146103	1
Caelo	Oak bark, cut	15146104	3
Caelo	Oak bark, cut	15146110	1
Caelo	Oak bark, cut	15146109	2
Caelo	Oak bark, cut	15146102	1
Caelo	Oak bark, cut	14277004	1
Caelo	Oak bark, cut	14277004	2
Caelo	Oak bark, cut	15146101	2
Caelo	Oak bark, cut	13360308	1
Caelo	Oak bark, cut	14277004	1
Caelo	Oak bark, cut	10061504	1
Caelo	Oak bark, cut	13360309	1
Caelo	Oak bark, cut	13360311	1
Caelo	Oak bark, cut	13360306	1
Caelo	Oak bark, cut	13360312	5
Caelo	Oak bark, cut	13360304	1
Caelo	Oak bark, cut	280814A	1
Caelo	Oak bark, cut	15314904	1
Caelo	Oak bark, cut	13360305	1
Caelo	Oak bark, cut	1081401	1
Caelo	Oak bark, cut	13360306	2
Caelo	Oak bark, cut	13360303	1
Caelo	Oak bark, cut	13360301	1
Caelo	Oak bark, cut	13360302	6
Caelo	Oak bark, cut	11351210	1
Caelo	Oak bark, cut	41109-4-P	1
Caelo	Oak bark, cut	11351210	8
Caelo	Oak bark, cut	15314911	1
Caelo	Oak bark, cut	16081305	1
Caelo	Oak bark, cut	12221105	2
Caelo	Oak bark, cut	11351201	3
Caelo	Oak bark, cut	19032013	2
Caelo	Oak bark, cut	15314905	2
Caelo	Oak bark, cut	15314903	1
Caelo	Oak bark, cut	15314910	2
Caelo	Oak bark, cut	15314901	1
Caelo	Oak bark, cut	9357419	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Oak bark, cut	12221104	2
Caelo	Oak bark, cut	10238612	1
Caelo	Oak bark, cut	10348612	1
Caelo	Oak bark, cut	11351202	5
Caelo	Oak bark, cut	7021302	1
Caelo	Oak bark, cut	12221101	9
Caelo	Oak bark, cut	11351207	11
Caelo	Oak bark, cut	5091302	1
Caelo	Oak bark, cut	10348615	1
Caelo	Oak bark, cut	26081302	1
Caelo	Oak bark, cut	11351205	3
Caelo	Oak bark, cut	12221102	1
Caelo	Oak bark, cut	15314908	1
Caelo	Oak bark, cut	3041301	1
Caelo	Oak bark, cut	22081204	1
Caelo	Oak bark, powdered	10261513	1
Caelo	Oat straw, cut	11280015	3
Caelo	Oat straw, cut	11280008	2
Caelo	Oat straw, cut	11280012	2
Caelo	Oat straw, cut	11280017	2
Caelo	Oat straw, cut	11280014	1
Caelo	Octenidine dihydrochl...	13166808	2
Caelo	Octenidine dihydrochl...	15275504	1
Caelo	Octenidine dihydrochl...	6325305	1
Caelo	Octenidine dihydrochl...	16325410	2
Caelo	Octenidine dihydrochl...	17042602	7
Caelo	Octenidine dihydrochl...	17042305	3
Caelo	Octenidine dihydrochl...	17042303	1
Caelo	Octenidine dihydrochl...	17042304	12
Caelo	Octenidine dihydrochl...	15006203	1
Caelo	Octenidine dihydrochl...	16325408	2
Caelo	Octenidine dihydrochl...	15006202	3
Caelo	Octenidine dihydrochl...	13166810	1
Caelo	Octenidine dihydrochl...	17042903	16
Caelo	Octenidine dihydrochl...	17042302	8
Caelo	Octenidine dihydrochl...	16325406	6
Caelo	Octenidine dihydrochl...	15275502	6
Caelo	Octenidine dihydrochl...	16302910	1
Caelo	Octenidine dihydrochl...	15356601	6
Caelo	Octenidine dihydrochl...	15356603	3
Caelo	Octenidine dihydrochl...	15356703	9
Caelo	Octenidine dihydrochl...	1515E-03732	1
Caelo	Octenidine dihydrochl...	16032803	7
Caelo	Octenidine dihydrochl...	16032901	2
Caelo	Octenidine dihydrochl...	15275503	2
Caelo	Octenidine dihydrochl...	816M-03732	1
Caelo	Octenidine dihydrochl...	16032902	4
Caelo	Octenidine dihydrochl...	152755020	1
Caelo	Octenidine dihydrochl...	16033002	2
Caelo	Octenidine dihydrochl...	16032909	7
Caelo	Octenidine dihydrochl...	16033001	2
Caelo	Octenidine dihydrochl...	16032910	5
Caelo	Octenidine dihydrochl...	16277802	10
Caelo	Octenidine dihydrochl...	16325304	10
Caelo	Octenidine dihydrochl...	16032903	1
Caelo	Octenidine dihydrochl...	16335505	4
Caelo	Octenidine dihydrochl...	13166801	1
Caelo	Octenidine dihydrochl...	16335504	16
Caelo	Octenidine dihydrochl...	1728506	1
Caelo	Octenidine dihydrochl...	16335506	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Octenidine dihydrochl...	16335502	3
Caelo	Octenidine dihydrochl...	18223419	1
Caelo	Octenidine dihydrochl...	17218051	1
Caelo	Octenidine dihydrochl...	17218507	8
Caelo	Octenidine dihydrochl...	11071925	1
Caelo	Octenidine dihydrochl...	18223407	4
Caelo	Octenidine dihydrochl...	18017402	33
Caelo	Octenidine dihydrochl...	17218503	8
Caelo	Octenidine dihydrochl...	13061801	1
Caelo	Octenidine dihydrochl...	300718A	1
Caelo	Octenidine dihydrochl...	17218907	1
Caelo	Octenidine dihydrochl...	18223403	16
Caelo	Octenidine dihydrochl...	18197501	18
Caelo	Octenidine dihydrochl...	18197502	9
Caelo	Octenidine dihydrochl...	18223418	1
Caelo	Octenidine dihydrochl...	18197504	7
Caelo	Octenidine dihydrochl...	2(18223413)	1
Caelo	Octenidine dihydrochl...	18223413	6
Caelo	Octenidine dihydrochl...	17218509	15
Caelo	Octenidine dihydrochl...	181975025789	1
Caelo	Octenidine dihydrochl...	17218508	18
Caelo	Octenidine dihydrochl...	172185035514	1
Caelo	Octenidine dihydrochl...	16335509	8
Caelo	Octenidine dihydrochl...	16335508	4
Caelo	Octenidine dihydrochl...	1711160HC/16335509	1
Caelo	Octenidine dihydrochl...	16335510	4
Caelo	Octenidine dihydrochl...	16335507	9
Caelo	Octenidine dihydrochl...	18223409	10
Caelo	Octenidine dihydrochl...	172118501	1
Caelo	Octenidine dihydrochl...	18223411	5
Caelo	Octenidine dihydrochl...	17218501	14
Caelo	Octenidine dihydrochl...	17218506	18
Caelo	Octenidine dihydrochl...	17042501	1
Caelo	Octenidine dihydrochl...	18223404	6
Caelo	Ofloxacin	15H19-B01	1
Caelo	Olive leaf	12159805	1
Caelo	Omeprazole	18E22-B03-353496	1
Caelo	Omeprazole	18K01-B03-358185	1
Caelo	Omeprazole	16120-B03-329040	1
Caelo	Oregano herb	12360901	6
Caelo	Oregano herb	11340709	1
Caelo	Oregano herb	11340703	2
Caelo	Oregano herb	12360902	3
Caelo	Oregano herb	11340708	1
Caelo	Oregano herb	12360904	1
Caelo	Oregano herb	11340705	1
Caelo	Orris root	13113402	1
Caelo	Orris root	11026706	1
Caelo	Orris root	13113401	2
Caelo	Orris root	1026706	1
Caelo	Orris root	11026704	2
Caelo	Orris root	13113404	1
Caelo	Oxalic acid	12212407	1
Caelo	Oxalic acid	7111808	1
Caelo	Oxalic acid	14334703	1
Caelo	Oxalic acid	13353915	1
Caelo	Oxalic acid	4812U-01039	1
Caelo	Oxalic acid	16050205	1
Caelo	Oxalic acid	17114003	1
Caelo	Oxalic acid	11111703	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Oxalic acid	16297401	1
Caelo	Oxalic acid	13195414	1
Caelo	Oxalic acid	13353901	1
Caelo	Oxalic acid	12147010	1
Caelo	Oxalic acid	13353907	2
Caelo	Oxalic acid	13353905	2
Caelo	Oxalic acid	12212411	1
Caelo	Oxalic acid	12212406	1
Caelo	Oxalic acid	13195413	1
Caelo	Oxalic acid	13353910	1
Caelo	Oxalic acid	10383104	1
Caelo	Oxalic acid	914E-01039	1
Caelo	Oxetacaine	17249810	2
Caelo	Oxybuprocaine hydroch...	14G21B02300384	1
Caelo	Oxybutynin hydrochlor...	15F05_B03_315221	5
Caelo	Oxybutynin hydrochlor...	16/23-B03-333371	3
Caelo	Oxybutynin hydrochlor...	15F05-B03-318524	4
Caelo	Oxybutynin hydrochlor...	16A18_B02_321181	2
Caelo	Oxybutynin hydrochlor...	18A12-B03-344350	1
Caelo	Oxybutynin hydrochlor...	16/23-B03-328916	1
Caelo	Oxybutynin hydrochlor...	16A18-B02-321181	3
Caelo	Oxybutynin hydrochlor...	18A12-B03	3
Caelo	Oxybutynin hydrochlor...	18A12-Bo3-344350	1
Caelo	Oxybutynin hydrochlor...	129099	2
Caelo	Oxytetracycline hydro...	1701034-02	27
Caelo	Oxytetracycline hydro...	17316812	2
Caelo	Oxytetracycline hydro...	14I22-B01-301463	1
Caelo	Oxytetracycline hydro...	16077012	1
Caelo	Oxytetracycline hydro...	16077017	1
Caelo	Oxytetracycline hydro...	1607716	1
Caelo	Oxytetracycline hydro...	163406010	1
Caelo	Oxytetracycline hydro...	16077015	1
Caelo	Oxytetracycline hydro...	16340605	4
Caelo	Oxytetracycline hydro...	17316806	4
Caelo	Oxytetracycline hydro...	17316805	5
Caelo	Oxytetracycline hydro...	16340609	1
Caelo	Oxytetracycline hydro...	16340610	12
Caelo	Oxytetracycline hydro...	16340606	3
Caelo	Oxytetracycline hydro...	16340611	2
Caelo	Oxytetracycline hydro...	1/01034-02	1
Caelo	Oxytetracycline hydro...	16077009	1
Caelo	Oxytetracycline hydro...	16077016	2
Caelo	Oxytetracycline hydro...	21021305	1
Caelo	Oxytetracycline hydro...	11336510	3
Caelo	Oxytetracycline hydro...	12328104	1
Caelo	Oxytetracycline hydro...	12328103	2
Caelo	Oxytetracycline hydro...	5041402	1
Caelo	Oxytetracycline hydro...	12328107	2
Caelo	Oxytetracycline hydro...	13L19-B33-290730	4
Caelo	Oxytetracycline hydro...	12091405	1
Caelo	Oxytetracycline hydro...	14128003	1
Caelo	Oxytetracycline hydro...	160202502	4
Caelo	Oxytetracycline hydro...	14128004	1
Caelo	Oxytetracycline hydro...	1407012-01	17
Caelo	Oxytetracycline hydro...	15060108	1
Caelo	Oxytetracycline hydro...	12130001	1
Caelo	Oxytetracycline hydro...	16077011	1
Caelo	Oxytetracycline hydro...	16077007	2
Caelo	Oxytetracycline hydro...	16044007	1
Caelo	Oxytetracycline hydro...	16077002	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Oxytetracycline hydro...	16077003	3
Caelo	Oxytetracycline hydro...	16077008	7
Caelo	Oxytetracycline hydro...	150603702	1
Caelo	Oxytetracycline hydro...	15060109	9
Caelo	Oxytetracycline hydro...	16077001	3
Caelo	Oxytetracycline hydro...	11336507	1
Caelo	Papaverine hydrochlor...	12316105	1
Caelo	Papaverine hydrochlor...	16315907	1
Caelo	Papaverine hydrochlor...	12166112	1
Caelo	Papaverine hydrochlor...	B:09157221	1
Caelo	Papaverine hydrochlor...	12166105	1
Caelo	Papaverine hydrochlor...	13324307	1
Caelo	Papaverine hydrochlor...	15160608	2
Caelo	Papaverine hydrochlor...	16315903	1
Caelo	Paracetamol	15217708	1
Caelo	Paracetamol	1814U-02045	1
Caelo	Paracetamol	15075704	4
Caelo	Paracetamol	13082010	2
Caelo	Paracetamol	13082011	1
Caelo	Paracetamol	13082005	1
Caelo	Paracetamol	15245902	1
Caelo	Paracetamol	10325020	1
Caelo	Paracetamol	1912M-02045	1
Caelo	Paracetamol	15217704	1
Caelo	Paracetamol	15075708	1
Caelo	Paracetamol	11300120	2
Caelo	Paracetamol	10065102	1
Caelo	Paracetamol	2815I-02045	1
Caelo	Paracetamol	15217710	3
Caelo	Paracetamol	15217712	1
Caelo	Paracetamol	18083005	1
Caelo	Paracetamol	18083003	1
Caelo	Paracetamol	17028411	1
Caelo	Paracetamol	11300108	1
Caelo	Paraffin, hard	13109203	4
Caelo	Paraffin, hard	14126202	1
Caelo	Paraffin, hard	1136206	1
Caelo	Paraffin, hard	11363220	3
Caelo	Paraffin, hard	13390403	1
Caelo	Paraffin, hard	13390413	1
Caelo	Paraffin, hard	13390406	4
Caelo	Paraffin, hard	13109224	5
Caelo	Paraffin, hard	12252205	2
Caelo	Paraffin, hard	13109221	1
Caelo	Paraffin, hard	11021202	1
Caelo	Paraffin, hard	11363205	1
Caelo	Paraffin, hard	13109214	4
Caelo	Paraffin, hard	12252206	1
Caelo	Paraffin, hard	14286503	1
Caelo	Paraffin, hard	13109213	6
Caelo	Paraffin, hard	12252211	3
Caelo	Paraffin, hard	10194101	2
Caelo	Paraffin, hard	12252221	3
Caelo	Paraffin, hard	13109218	1
Caelo	Paraffin, hard	14286504	1
Caelo	Paraffin, hard	123654	1
Caelo	Paraffin, hard	13390415	3
Caelo	Paraffin, hard	18070107	7
Caelo	Paraffin, hard	18070112	1
Caelo	Paraffin, hard	18070118	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Paraffin, hard	17245602	6
Caelo	Paraffin, hard	16277710	2
Caelo	Paraffin, hard	16277706	3
Caelo	Paraffin, hard	16274203	1
Caelo	Paraffin, hard	16277704	3
Caelo	Paraffin, hard	16000921	3
Caelo	Paraffin, hard	16000922	3
Caelo	Paraffin, hard	16000912	4
Caelo	Paraffin, hard	16000911	4
Caelo	Paraffin, hard	16000910	1
Caelo	Paraffin, hard	14286502	1
Caelo	Paraffin, hard	302709	1
Caelo	Paraffin, hard	15295407	1
Caelo	Paraffin, hard	15295414	3
Caelo	Paraffin, hard	15295410	5
Caelo	Paraffin, hard	15295406	2
Caelo	Paraffin, hard	15232001	1
Caelo	Paraffin, hard	15295405	5
Caelo	Paraffin, hard	152320004	1
Caelo	Paraffin, hard	15232004	5
Caelo	Paraffin, hard	1523002	1
Caelo	Paraffin, hard	13390408	1
Caelo	Paraffin, hard	17245604	1
Caelo	Paraffin, hard	14286511	2
Caelo	Paraffin, hard	14286517	5
Caelo	Paraffin, hard	16000903	2
Caelo	Paraffin, hard	17107304	1
Caelo	Paraffin, hard	17347701	13
Caelo	Paraffin, microcrysta...	15412808	1
Caelo	Paraffin, microcrysta...	12355001	1
Caelo	Paraffin, microcrysta...	13274301	2
Caelo	Paraffin, microcrysta...	13274303	5
Caelo	Paraffin, microcrysta...	14307110	1
Caelo	Paraffin, microcrysta...	14307108	2
Caelo	Paraffin, microcrysta...	7111201	1
Caelo	Paraffin, microcrysta...	52622395	1
Caelo	Paraffin, microcrysta...	15412801	2
Caelo	Paraffin, microcrysta...	12095101	1
Caelo	Paraffin, microcrysta...	30101202	1
Caelo	Paraffin, microcrysta...	16000911	1
Caelo	Paraffin, microcrysta...	7111202	1
Caelo	Paraffin, microcrysta...	3042012M1	1
Caelo	Parsley herb, rubbed	12376601	2
Caelo	Parsley herb, rubbed	12376602	2
Caelo	Parsley herb, rubbed	51978126	1
Caelo	Parsley herb, rubbed	12376603	2
Caelo	Parsley root, cut	12080309	1
Caelo	Parsley root, cut	12080304	1
Caelo	Parsley root, cut	13156003	1
Caelo	Parsley root, cut	13156005	1
Caelo	Parsley seed	PrüfNr. 2814A-01550	1
Caelo	Passionflower herb, cut	15085105	1
Caelo	Passionflower herb, cut	13141205	4
Caelo	Passionflower herb, cut	13141203	2
Caelo	Passionflower herb, cut	11049111	1
Caelo	Passionflower herb, cut	13141201	1
Caelo	Passionflower herb, cut	12031712	1
Caelo	Passionflower herb, cut	11049108	1
Caelo	Passionflower herb, cut	12031716	3
Caelo	Passionflower herb, cut	12031704	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Passionflower herb, cut	12031706	3
Caelo	Passionflower herb, cut	13141202	5
Caelo	Passionflower herb, cut	12031714	2
Caelo	Passionflower herb, cut	12031701	1
Caelo	Passionflower herb, cut	12031715	7
Caelo	Passionflower herb, cut	12031717	9
Caelo	Passionflower herb, cut	12101301	1
Caelo	Passionflower herb, cut	1011A-01669	1
Caelo	Passionflower herb, cut	12031705	4
Caelo	Passionflower herb, cut	11049116	4
Caelo	Passionflower herb, cut	12031703	2
Caelo	Passionflower herb, cut	12031710	13
Caelo	Patent blue V	13351803	1
Caelo	Patent blue V	17034204	4
Caelo	Patent blue V	12364202	1
Caelo	Patent blue V	16306904	1
Caelo	Patent blue V	15060204	1
Caelo	Patent blue V	12364206	1
Caelo	Patent blue V	12364206	1
Caelo	Patent blue V	17034207	1
Caelo	Patent blue V	11305303	3
Caelo	Patent blue V	150602	1
Caelo	Patent blue V	13351807	5
Caelo	Patent blue V	13351808	1
Caelo	Patent blue V	11305302	2
Caelo	Peppermint leaf, cut	13030804	1
Caelo	Peppermint leaf, cut	161212	1
Caelo	Peppermint leaf, cut	Int. Ch. : 313011	1
Caelo	Peppermint leaf, cut	ChB: 311092	1
Caelo	Peppermint leaf, cut	12054149	1
Caelo	Peppermint leaf, cut	12266903	5
Caelo	Peppermint leaf, cut	11347609	1
Caelo	Peppermint leaf, cut	13105001	1
Caelo	Peppermint leaf, cut	13030806	1
Caelo	Peppermint leaf, cut	12266904	2
Caelo	Peppermint leaf, cut	13030808	1
Caelo	Peppermint leaf, cut	13030808-3401310161	1
Caelo	Peppermint leaf, cut	5091307	1
Caelo	Peppermint leaf, cut	12266913	9
Caelo	Peppermint leaf, cut	2181B090928	1
Caelo	Peppermint leaf, cut	8101301	1
Caelo	Peppermint leaf, cut	2181A1209-13	1
Caelo	Peppermint leaf, cut	13030809	2
Caelo	Peppermint leaf, cut	13030807	3
Caelo	Peppermint leaf, cut	12032311	1
Caelo	Peppermint leaf, cut	1132803	1
Caelo	Peppermint leaf, cut	13214720	1
Caelo	Peppermint leaf, cut	12266912	1
Caelo	Peppermint leaf, cut	13030805	5
Caelo	Peppermint leaf, cut	11294401	3
Caelo	Peppermint leaf, cut	11294409	2
Caelo	Peppermint leaf, cut	268320	1
Caelo	Peppermint leaf, cut	11327801	2
Caelo	Peppermint leaf, cut	20061202	1
Caelo	Peppermint leaf, cut	11327803	5
Caelo	Peppermint leaf, cut	13214705	1
Caelo	Peppermint leaf, cut	12032319	2
Caelo	Peppermint leaf, cut	20111215	1
Caelo	Peppermint leaf, cut	14083102	2
Caelo	Peppermint leaf, cut	14339202	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Peppermint leaf, cut	13214707	2
Caelo	Peppermint leaf, cut	14221701	1
Caelo	Peppermint leaf, cut	6011403	1
Caelo	Peppermint leaf, cut	13214706	1
Caelo	Peppermint leaf, cut	12266908	7
Caelo	Phenazone	18061101	1
Caelo	Pheniramine maleate	14C25-B01-302272	1
Caelo	Pheniramine maleate	14C25-B01-322989	1
Caelo	Pheniramine maleate	12352604	1
Caelo	Pheniramine maleate	15172704	1
Caelo	Pheniramine maleate	9276210	1
Caelo	Pheniramine maleate	15172711	1
Caelo	Pheniramine maleate	12352602	1
Caelo	Pheniramine maleate	17316702	1
Caelo	Pheniramine maleate	15172708	1
Caelo	Pheniramine maleate	1235611	1
Caelo	Phenol	11163802	1
Caelo	Phenol	12223020	1
Caelo	Phenol	12223014	1
Caelo	Phenol	7061305	1
Caelo	Phenol	11163809	1
Caelo	Phenol	11163801	1
Caelo	Phenol	12223001	1
Caelo	Phenol	12223022	1
Caelo	Phenylbutazon	10262901	1
Caelo	Phenylbutazon	15218502	1
Caelo	Phenylephrine hydroch...	14H13-B02-308438	1
Caelo	Phenytoin	17L22-B04-345861	1
Caelo	Pilocarpine hydrochlo...	12K07-N09	1
Caelo	Pilocarpine hydrochlo...	18138704	2
Caelo	Pilocarpine hydrochlo...	2211112032607	1
Caelo	Pilocarpine hydrochlo...	17236401	1
Caelo	Pilocarpine hydrochlo...	17124208	3
Caelo	Pilocarpine hydrochlo...	12217316	1
Caelo	Pilocarpine hydrochlo...	12217307	1
Caelo	Pilocarpine hydrochlo...	25051303	1
Caelo	Pilocarpine hydrochlo...	12217302	2
Caelo	Piper cubeba	11081903	1
Caelo	Podophyllum resin (fr...	11085807	4
Caelo	Podophyllum resin (fr...	12287810	1
Caelo	Podophyllum resin (fr...	12379903	1
Caelo	Podophyllum resin (fr...	82433349	1
Caelo	Podophyllum resin (fr...	11085811	1
Caelo	Podophyllum resin (fr...	12287805	1
Caelo	Poloxamers 407	13133201	1
Caelo	Poloxamers 407	13133203	2
Caelo	Poloxamers 407	12286806	1
Caelo	Poloxamers 407	18160702	3
Caelo	Poloxamers 407	17048601	3
Caelo	Poloxamers 407	12286805	1
Caelo	Poloxamers 407	13133207	1
Caelo	Poloxamers 407	12286801	2
Caelo	Poloxamers 407	18025101	7
Caelo	Poloxamers 407	18160704	1
Caelo	Poloxamers 407	17048606	2
Caelo	Poloxamers 407	18263605	1
Caelo	Poloxamers 407	16242605	1
Caelo	Poloxamers 407	14196101	4
Caelo	Poloxamers 407	162426056	1
Caelo	Poloxamers 407	17292101	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Poloxamers 407	14196108	2
Caelo	Poloxamers 407	14196105	1
Caelo	Poloxamers 407	14359501	1
Caelo	Poloxamers 407	16189301	1
Caelo	Poloxamers 407	16242601	3
Caelo	Poloxamers 407	16189302	4
Caelo	Poloxamers 407	15397802	6
Caelo	Poloxamers 407	14196104	1
Caelo	Polymyxin B sulfate	15350101	5
Caelo	Polymyxin B sulfate	13H28-B02-287953	1
Caelo	Polymyxin B sulfate	1102013-01	1
Caelo	Polymyxin B sulfate	1604016-01	3
Caelo	Polymyxin B sulfate	150725	1
Caelo	Polymyxin B sulfate	1802006-011	1
Caelo	Potassium citrate	15158706	1
Caelo	Potassium citrate	15148707	4
Caelo	Potassium citrate	13135907	4
Caelo	Potassium citrate	11320201	1
Caelo	Potassium citrate	13135902	1
Caelo	Potassium citrate	5555	1
Caelo	Potassium citrate	3112M01737	1
Caelo	Potassium citrate		12
Caelo	Potassium citrate	12345678	1
Caelo	Potassium citrate	11320212	3
Caelo	Potassium citrate	1112008-01	1
Caelo	Potassium citrate	123654879	1
Caelo	Potassium citrate	12135302	6
Caelo	Potassium citrate	28021402	1
Caelo	Potassium citrate	3456	1
Caelo	Potassium citrate	1234	6
Caelo	Potassium citrate	15123103	5
Caelo	Potassium citrate	13317602	2
Caelo	Potassium citrate	14138118	2
Caelo	Potassium citrate	189645	1
Caelo	Potassium citrate	5021509	1
Caelo	Potassium citrate	14138107	1
Caelo	Potassium citrate	15123102	2
Caelo	Potassium citrate	13317618	2
Caelo	Potassium citrate	15148701	7
Caelo	Potassium citrate	13317604	2
Caelo	Potassium citrate	15148706	9
Caelo	Potassium citrate	123test	2
Caelo	Potassium citrate	17061402	1
Caelo	Potassium citrate	2114M1737	1
Caelo	Potassium citrate	15148702	1
Caelo	Potassium citrate	13317613	3
Caelo	Potassium citrate	123Test	3
Caelo	Potassium citrate	12334	1
Caelo	Potassium citrate	123Testmessung	1
Caelo	Potassium citrate	13135908	1
Caelo	Potassium citrate	13317603	2
Caelo	Potassium citrate	15435701	6
Caelo	Potassium citrate	71730367	1
Caelo	Potassium citrate	154357	3
Caelo	Potassium citrate	16113003	1
Caelo	Potassium citrate	232262310107	1
Caelo	Potassium citrate	16113002	2
Caelo	Potassium citrate	16113006	4
Caelo	Potassium citrate	13135904	1
Caelo	Potassium citrate	15435702	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Potassium citrate	16113012	3
Caelo	Potassium citrate	16113008	1
Caelo	Potassium citrate	16113007	1
Caelo	Potassium citrate	13080603	1
Caelo	Potassium citrate	17092904	2
Caelo	Potassium citrate	18007603	1
Caelo	Potassium citrate	18038506	3
Caelo	Potassium citrate	18038502	1
Caelo	Potassium citrate	1543702	1
Caelo	Potassium citrate	gfte	1
Caelo	Potassium citrate	15435710	5
Caelo	Potassium citrate	3012Q-01737	1
Caelo	Potassium citrate	12135303	2
Caelo	Potassium citrate	12648	1
Caelo	Potassium citrate	9259107	1
Caelo	Potassium citrate	12135309	2
Caelo	Potassium citrate	1215309	1
Caelo	Potassium citrate	12299201	5
Caelo	Potassium citrate	12299203	4
Caelo	Potassium citrate	11102201	1
Caelo	Potassium citrate	71730128	1
Caelo	Potassium citrate	12299202	1
Caelo	Potassium citrate	9259120	1
Caelo	Potassium citrate	411025-1-P	1
Caelo	Potassium citrate	12456	1
Caelo	Potassium citrate	14138105	6
Caelo	Potassium citrate	151231	1
Caelo	Potassium citrate	154357025	1
Caelo	Potassium citrate	13080602	8
Caelo	Potassium citrate	14138119	1
Caelo	Potassium hydrocarbon...	12141924	3
Caelo	Potassium hydrocarbon...	20080410-3	1
Caelo	Potassium hydrocarbon...	11105822	1
Caelo	Potassium hydrocarbon...	10052610	1
Caelo	Potassium hydrocarbon...	10052622	1
Caelo	Potassium hydrocarbon...	12141905	3
Caelo	Potassium hydrocarbon...	12239703	3
Caelo	Potassium hydrocarbon...	13390307	1
Caelo	Potassium hydrocarbon...	12141930	2
Caelo	Potassium hydrocarbon...	12141912	3
Caelo	Potassium hydrocarbon...	24051301	1
Caelo	Potassium hydrocarbon...	13390303	1
Caelo	Potassium hydrocarbon...	10052629	1
Caelo	Potassium hydrocarbon...	12141915	4
Caelo	Potassium hydrocarbon...	13141926	1
Caelo	Potassium hydrocarbon...	IntCh:217002	1
Caelo	Potassium hydrocarbon...	1303002-01	5
Caelo	Potassium hydrocarbon...	12141920	1
Caelo	Potassium hydrocarbon...	1214190	1
Caelo	Potassium hydrocarbon...	12141919	4
Caelo	Potassium hydrocarbon...	700544	1
Caelo	Potassium hydrocarbon...	3812Q-01728	1
Caelo	Potassium hydrocarbon...	1214924	1
Caelo	Potassium hydrocarbon...	11105811	3
Caelo	Potassium hydrocarbon...	1011M-01728	1
Caelo	Potassium hydrocarbon...	12141910	8
Caelo	Potassium hydrocarbon...	8101202	1
Caelo	Potassium hydrocarbon...	12239702	2
Caelo	Potassium hydrocarbon...	11105818	1
Caelo	Potassium hydrocarbon...	11105804	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Potassium hydrocarbon...	11105812	1
Caelo	Potassium hydrocarbon...	11105808	1
Caelo	Potassium hydrocarbon...	12239701	3
Caelo	Potassium hydrocarbon...	6011405	1
Caelo	Potassium hydrogen ta...	15354614	2
Caelo	Potassium hydrogen ta...	23062015A	1
Caelo	Potassium hydrogen ta...	14115208	4
Caelo	Potassium hydrogen ta...	15081008	2
Caelo	Potassium hydrogen ta...	90573259	1
Caelo	Potassium hydrogen ta...	13133104	2
Caelo	Potassium hydrogen ta...	15354613	2
Caelo	Potassium hydrogen ta...	11127417	1
Caelo	Potassium hydrogen ta...	1534609	1
Caelo	Potassium hydrogen ta...	13030006	2
Caelo	Potassium hydrogen ta...	13133114	1
Caelo	Potassium hydrogen ta...	15354601	1
Caelo	Potassium hydrogen ta...	15081010	1
Caelo	Potassium hydrogen ta...	17061303	1
Caelo	Potassium hydrogen ta...	11345807	1
Caelo	Potassium hydrogen ta...	13133120	2
Caelo	Potassium hydrogen ta...	160514b	1
Caelo	Potassium hydrogen ta...	15081005	1
Caelo	Potassium hydroxide	11165708	1
Caelo	Potassium hydroxide	17091202	1
Caelo	Potassium hydroxide	12057912	2
Caelo	Potassium hydroxide	12382611	1
Caelo	Potassium hydroxide	12382601	2
Caelo	Potassium hydroxide	12057908	2
Caelo	Potassium hydroxide	12382603	4
Caelo	Potassium hydroxide	27111310	1
Caelo	Potassium hydroxide	12057906	1
Caelo	Potassium hydroxide	13275207	1
Caelo	Potassium hydroxide	12382606	2
Caelo	Potassium hydroxide	13275201	2
Caelo	Potassium hydroxide	11165705	1
Caelo	Potassium hydroxide	12057905	1
Caelo	Potassium sodium tart...	15375601	12
Caelo	Potassium sodium tart...	17257701	3
Caelo	Potassium sorbate	12211108	2
Caelo	Potassium sorbate	7111805	1
Caelo	Potassium sorbate	18165006	6
Caelo	Potassium sorbate	14380601	19
Caelo	Potassium sorbate	12059908	2
Caelo	Potassium sorbate	18165010	9
Caelo	Potassium sorbate	13224610	1
Caelo	Potassium sorbate	1816510	1
Caelo	Potassium sorbate	12059909	4
Caelo	Potassium sorbate	14192605	2
Caelo	Potassium sorbate	14380606	8
Caelo	Potassium sorbate	1438060	2
Caelo	Potassium sorbate	AKD10	1
Caelo	Potassium sorbate	14380602	2
Caelo	Potassium sorbate	14380610	13
Caelo	Potassium sorbate	14380607	5
Caelo	Potassium sorbate	13224612	18
Caelo	Potassium sorbate	13111409	1
Caelo	Potassium sorbate	17345116	13
Caelo	Potassium sorbate	17170007	16
Caelo	Potassium sorbate	9305601	2
Caelo	Potassium sorbate	17345103	14

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Potassium sorbate	17345105	1
Caelo	Potassium sorbate	13065217	15
Caelo	Potassium sorbate	12059902	1
Caelo	Potassium sorbate	200613K1	1
Caelo	Potassium sorbate	12345	1
Caelo	Potassium sorbate	13065211	1
Caelo	Potassium sorbate	9061806	1
Caelo	Potassium sorbate	11192305	3
Caelo	Potassium sorbate	90318101	1
Caelo	Potassium sorbate	13277804	1
Caelo	Potassium sorbate	17170017	11
Caelo	Potassium sorbate	13224611	14
Caelo	Potassium sorbate	13224608	2
Caelo	Potassium sorbate	17345104	20
Caelo	Potassium sorbate	14380618	1
Caelo	Potassium sorbate	14380606,PZN01397235	1
Caelo	Potassium sorbate	16326602	2
Caelo	Potassium sorbate	12211121	4
Caelo	Potassium sorbate	14380627	10
Caelo	Potassium sorbate	17170013	18
Caelo	Potassium sorbate	15382105	11
Caelo	Potassium sorbate	12059914	2
Caelo	Potassium sorbate	16326615	1
Caelo	Potassium sorbate	14380620	6
Caelo	Potassium sorbate	14380632	20
Caelo	Potassium sorbate	17170007-20171106Ma1	1
Caelo	Potassium sorbate	13065205	13
Caelo	Potassium sorbate	12211102	4
Caelo	Potassium sorbate	17170005	1
Caelo	Potassium sorbate	13065204	12
Caelo	Potassium sorbate	17170002	4
Caelo	Potassium sorbate	K-010213-1	1
Caelo	Potassium sorbate	17170006	9
Caelo	Potassium sorbate	14140007	1
Caelo	Potassium sorbate	91623299	2
Caelo	Potassium sorbate	16326606	13
Caelo	Potassium sorbate	14380613	1
Caelo	Potassium sorbate	16326609	12
Caelo	Potassium sorbate	14380612	2
Caelo	Potassium sorbate	15248204	17
Caelo	Potassium sorbate	143806	1
Caelo	Potassium sorbate	17170012	2
Caelo	Potassium sorbate	15382106	1
Caelo	Potassium sorbate	15248205	11
Caelo	Potassium sorbate	14380622	12
Caelo	Potassium sorbate	12211107	16
Caelo	Potassium sorbate	14380621	10
Caelo	Potassium sorbate	16326608	16
Caelo	Potassium sorbate	15248203	3
Caelo	Potassium sorbate	1205999	1
Caelo	Potassium sorbate	12059905	2
Caelo	Potassium sorbate	12211101	10
Caelo	Potassium sorbate	3512M-01748	1
Caelo	Potassium sorbate	20021112-2	1
Caelo	Potassium sorbate	1112331024	1
Caelo	Potassium sorbate	14380631	1
Caelo	Potassium sorbate	12211123	1
Caelo	Potassium sorbate	16326607	18
Caelo	Potassium sorbate	16156505	1
Caelo	Potassium sorbate	14380633	14

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Potassium sorbate	14380635	7
Caelo	Potassium sorbate	15382114	10
Caelo	Potassium sorbate	11214207	1
Caelo	Potassium sorbate	27101705	1
Caelo	Potassium sorbate	10012210	1
Caelo	Potassium sorbate	15382103	1
Caelo	Potassium sorbate	10199505	1
Caelo	Potassium tartrate he...	19051804	1
Caelo	Potassium tartrate he...	11225007	1
Caelo	Potassium tartrate he...	12322110	1
Caelo	Potassium tartrate he...	13111415	1
Caelo	Potassium tartrate he...	13331113	1
Caelo	Potassium tartrate he...	15167703	1
Caelo	Potassium tartrate he...	27061304	1
Caelo	Potassium tartrate he...	17176002	2
Caelo	Potassium tartrate he...	17176004	1
Caelo	Potassium tartrate he...	17176007	1
Caelo	Potassium tartrate he...	24061705	1
Caelo	Potato starch	15062912	2
Caelo	Potato starch	14033815	1
Caelo	Potato starch	11301012	2
Caelo	Potato starch	17327201	3
Caelo	Povidone 25, soluble	14207517	1
Caelo	Povidone 25, soluble	17209103	1
Caelo	Povidone 25, soluble	14207515	1
Caelo	Povidone 25, soluble	18024602	1
Caelo	Povidone 25, soluble	18027602	2
Caelo	Povidone 25, soluble	15383103	1
Caelo	Povidone 25, soluble	18112201	4
Caelo	Povidone 25, soluble	16051102	1
Caelo	Povidone 25, soluble	17209104	3
Caelo	Povidone 25, soluble	14207511	1
Caelo	Povidone 25, soluble	15383102	1
Caelo	Povidone 25, soluble	16051104	1
Caelo	Povidone 25, soluble	16051101	2
Caelo	Povidone 25, soluble	16307812	1
Caelo	Povidone 25, soluble	16307809	3
Caelo	Povidone 25, soluble	14142103	1
Caelo	Povidone 25, soluble	14207508	3
Caelo	Povidone 25, soluble	16051108	1
Caelo	Povidone 25, soluble	163078085328	1
Caelo	Povidone 25, soluble	16051110	1
Caelo	Povidone 25, soluble	14207511	3
Caelo	Povidone 25, soluble	14142109	1
Caelo	Povidone 25, soluble	16307811	1
Caelo	Povidone 25, soluble	18182002	2
Caelo	Povidone 25, soluble	15383101	3
Caelo	Povidone 25, soluble	12032714	1
Caelo	Povidone 25, soluble	14142109	2
Caelo	Povidone 25, soluble	12376301	1
Caelo	Povidone 25, soluble	13323407	2
Caelo	Povidone iodine	12153119	2
Caelo	Povidone iodine	12153110	2
Caelo	Povidone iodine	12153116	2
Caelo	Povidone iodine	19061302	1
Caelo	Povidone iodine	12072605	1
Caelo	Povidone iodine	17028102	2
Caelo	Povidone iodine	11168607	1
Caelo	Povidone iodine	15234109	1
Caelo	Povidone iodine	12153106	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Povidone iodine	15442405	1
Caelo	Povidone iodine	15234101	1
Caelo	Povidone iodine	14195810	2
Caelo	Povidone iodine	15121402	1
Caelo	Prasterone	16102-B02-326485	2
Caelo	Prasterone	16102-B02-326485	1
Caelo	Prasterone	16B19-B04-318665	4
Caelo	Prednicarbate, micron...	17220407	23
Caelo	Prednicarbate, micron...	C0333/1111070	1
Caelo	Prednicarbate, micron...	17220406	2
Caelo	Prednicarbate, micron...	17220408	25
Caelo	Prednicarbate, micron...	17014006	1
Caelo	Prednicarbate, micron...	15446302	10
Caelo	Prednicarbate, micron...	17315603	2
Caelo	Prednicarbate, micron...	172204075515	1
Caelo	Prednicarbate, micron...	15446305	6
Caelo	Prednicarbate, micron...	17220401	3
Caelo	Prednicarbate, micron...	17014003	15
Caelo	Prednicarbate, micron...	160500	1
Caelo	Prednicarbate, micron...	1517E-03595	1
Caelo	Prednicarbate, micron...	17013804	15
Caelo	Prednicarbate, micron...	17013805	6
Caelo	Prednicarbate, micron...	16087411	3
Caelo	Prednicarbate, micron...	16087412	5
Caelo	Prednicarbate, micron...	16A18-B01-322733	2
Caelo	Prednicarbate, micron...	17014001	2
Caelo	Prednicarbate, micron...	17014005	12
Caelo	Prednicarbate, micron...	1608710	1
Caelo	Prednicarbate, micron...	16987410	1
Caelo	Prednicarbate, micron...	16087410	14
Caelo	Prednicarbate, micron...	1704003	1
Caelo	Prednicarbate, micron...	16A18-B01-322732	4
Caelo	Prednicarbate, micron...	15446304	1
Caelo	Prednicarbate, micron...	15446303	2
Caelo	Prednicarbate, micron...	16087408	20
Caelo	Prednicarbate, micron...	1704005	1
Caelo	Prednicarbate, micron...	15J22-B-321005	1
Caelo	Prednicarbate, micron...	16087401	13
Caelo	Prednicarbate, micron...	15446306	7
Caelo	Prednicarbate, micron...	2(18194104)	1
Caelo	Prednicarbate, micron...	18197115	1
Caelo	Prednicarbate, micron...	1507011-01	3
Caelo	Prednicarbate, micron...	14J16-B06	3
Caelo	Prednicarbate, micron...	15D15-B03	1
Caelo	Prednicarbate, micron...	15181403	11
Caelo	Prednicarbate, micron...	13224401	6
Caelo	Prednicarbate, micron...	12383309	1
Caelo	Prednicarbate, micron...	18194115	12
Caelo	Prednicarbate, micron...	18194112	2
Caelo	Prednicarbate, micron...	13372419	4
Caelo	Prednicarbate, micron...	17315607	2
Caelo	Prednicarbate, micron...	18194103	3
Caelo	Prednicarbate, micron...	9764897	1
Caelo	Prednicarbate, micron...	13372410	8
Caelo	Prednicarbate, micron...	13119401	1
Caelo	Prednicarbate, micron...	17315609	14
Caelo	Prednicarbate, micron...	13224406	1
Caelo	Prednicarbate, micron...	13119402	1
Caelo	Prednicarbate, micron...	18194105	22
Caelo	Prednicarbate, micron...	18194117	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Prednicarbate, micron...	13372403	10
Caelo	Prednicarbate, micron...	17220410	1
Caelo	Prednicarbate, micron...	19051601	8
Caelo	Prednicarbate, micron...	3(19051602)	1
Caelo	Prednicarbate, micron...	19051603	1
Caelo	Prednicarbate, micron...	15181409	1
Caelo	Prednicarbate, micron...	12383301	1
Caelo	Prednicarbate, micron...	Ch.B.:17220408	1
Caelo	Prednicarbate, micron...	18194104	14
Caelo	Prednicarbate, micron...	15181404	2
Caelo	Prednicarbate, micron...	15446301	13
Caelo	Prednicarbate, micron...	15181411	3
Caelo	Prednicarbate, micron...	18B27-.B01	1
Caelo	Prednicarbate, micron...	17315605	10
Caelo	Prednicarbate, micron...	19051604	1
Caelo	Prednisolon dihydroge...	14196602	1
Caelo	Prednisolon dihydroge...	17077502	2
Caelo	Prednisolon dihydroge...	18194011	1
Caelo	Prednisolon dihydroge...	16174702	1
Caelo	Prednisolon dihydroge...	18194006	2
Caelo	Prednisolon dihydroge...	18194001	1
Caelo	Prednisolon dihydroge...	18024604	1
Caelo	Prednisolon dihydroge...	16174701	1
Caelo	Prednisolon dihydroge...	16174710	1
Caelo	Prednisolon dihydroge...	16174706	1
Caelo	Prednisolon dihydroge...	15199602	2
Caelo	Prednisolon dihydroge...	18024602	1
Caelo	Prednisolon dihydroge...	15199603	1
Caelo	Prednisolone acetate	27051301	1
Caelo	Prednisolone acetate	17185603	14
Caelo	Prednisolone acetate	12091807	6
Caelo	Prednisolone acetate	12091804	8
Caelo	Prednisolone acetate	12091803	3
Caelo	Prednisolone acetate	14175513	13
Caelo	Prednisolone acetate	12324503	1
Caelo	Prednisolone acetate	14175508	14
Caelo	Prednisolone acetate	14175514	7
Caelo	Prednisolone acetate	1480905	1
Caelo	Prednisolone acetate	18067302	6
Caelo	Prednisolone acetate	14175509	9
Caelo	Prednisolone acetate	15055201	9
Caelo	Prednisolone acetate	15055206	9
Caelo	Prednisolone acetate	15055203	2
Caelo	Prednisolone acetate	18004402	6
Caelo	Prednisolone acetate	14380901	15
Caelo	Prednisolone acetate	14380904	2
Caelo	Prednisolone acetate	15055202	18
Caelo	Prednisolone acetate	14175512	1
Caelo	Prednisolone acetate	17086101	9
Caelo	Prednisolone acetate	12091802	8
Caelo	Prednisolone acetate	18067301	2
Caelo	Prednisolone acetate	11158003	1
Caelo	Prednisolone acetate	17185604	2
Caelo	Prednisolone acetate	11158002	1
Caelo	Prednisolone acetate	14380909	13
Caelo	Prednisolone acetate	17086106	2
Caelo	Prednisolone acetate	12091808	5
Caelo	Prednisolone acetate	1406014-02	15
Caelo	Prednisolone acetate	14380302	1
Caelo	Prednisolone acetate	17086105	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Prednisolone acetate	13262603	8
Caelo	Prednisolone acetate	13262605	14
Caelo	Prednisolone acetate	12324509	2
Caelo	Prednisolone acetate	14380905	16
Caelo	Prednisolone acetate	1234506	1
Caelo	Prednisolone acetate	4313A-02106	1
Caelo	Prednisolone acetate	13262606	4
Caelo	Prednisolone acetate	17185602	14
Caelo	Prednisolone acetate	12324508	9
Caelo	Prednisolone acetate	12324506	13
Caelo	Prednisolone acetate	11080403	1
Caelo	Prednisolone acetate	12324505	16
Caelo	Prednisolone acetate	12324514	5
Caelo	Prednisolone acetate	12324504	6
Caelo	Prednisolone acetate		2
Caelo	Prednisolone acetate	313I-02106	1
Caelo	Prednisolone acetate	12091810	2
Caelo	Prednisolone acetate	12324502	12
Caelo	Prednisolone acetate	12091815	2
Caelo	Prednisolone acetate	17086104	13
Caelo	Prednisolone acetate	18004903	1
Caelo	Prednisolone acetate	12091813	1
Caelo	Prednisolone acetate	12091806	3
Caelo	Prednisolone acetate	12091814	7
Caelo	Prednisolone acetate	11L21-N03	1
Caelo	Prednisolone acetate	6021306	1
Caelo	Prednisolone acetate	13262602	8
Caelo	Prednisolone acetate	13262601	15
Caelo	Prednisolone acetate	162468	2
Caelo	Prednisolone acetate	14380910	13
Caelo	Prednisolone acetate	15398402	12
Caelo	Prednisolone acetate	170104101	4
Caelo	Prednisolone acetate	15398407	8
Caelo	Prednisolone acetate	170861015494	1
Caelo	Prednisolone acetate	162446808	1
Caelo	Prednisolone acetate	17086103/1110917	1
Caelo	Prednisolone acetate	16246809	1
Caelo	Prednisolone acetate	18398303	1
Caelo	Prednisolone acetate	32304015	1
Caelo	Prednisolone acetate	1867314	1
Caelo	Prednisolone acetate	18067307	8
Caelo	Prednisolone acetate	18067317	5
Caelo	Prednisolone acetate	1624806	1
Caelo	Prednisolone acetate	15344002	1
Caelo	Prednisolone acetate	32304	1
Caelo	Prednisolone acetate	16246801	14
Caelo	Prednisolone acetate	16246804	10
Caelo	Prednisolone acetate	16246802	9
Caelo	Prednisolone acetate	15396304	1
Caelo	Prednisolone acetate	16246803	20
Caelo	Prednisolone acetate	160701001	4
Caelo	Prednisolone acetate	14380903	2
Caelo	Prednisolone acetate	717Q-02106	1
Caelo	Prednisolone acetate	165246802	1
Caelo	Prednisolone acetate	16246808	11
Caelo	Prednisolone acetate	18067309-P20180115Ma	1
Caelo	Prednisolone acetate	16246807	14
Caelo	Prednisolone acetate	18067310	3
Caelo	Prednisolone acetate	180500903	1
Caelo	Prednisolone acetate	16246805	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Prednisolone acetate	1624801	1
Caelo	Prednisolone acetate	16246811	4
Caelo	Prednisolone acetate	14380902	12
Caelo	Prednisolone acetate	15398307	9
Caelo	Prednisolone acetate	16248808	1
Caelo	Prednisolone acetate	180673025724	1
Caelo	Prednisolone acetate	17086102	7
Caelo	Prednisolone acetate	1616I-02106	1
Caelo	Prednisolone acetate	15398302	13
Caelo	Prednisolone acetate	18067309	12
Caelo	Prednisolone acetate	15398304	6
Caelo	Prednisolone acetate	18004401	11
Caelo	Prednisolone acetate	15398303	14
Caelo	Prednisolone acetate	15398401	9
Caelo	Prednisolone acetate	17086103	16
Caelo	Prednisolone acetate	1475509	1
Caelo	Prednisolone acetate	171856010	1
Caelo	Prednisolone acetate	15398306	6
Caelo	Prednisolone acetate	20160727Ma1	1
Caelo	Prednisolone hydrogen...	15C11-B05-305533	1
Caelo	Prednisolone hydrogen...	12A20-N11	1
Caelo	Prednisolone hydrogen...	19081502	1
Caelo	Prednisolone hydrogen...	15C11B05328275	1
Caelo	Prednisolone, microni...	13343104	3
Caelo	Prednisolone, microni...	152006105	1
Caelo	Prednisolone, microni...	1234556	1
Caelo	Prednisolone, microni...	1212021-03	2
Caelo	Prednisolone, microni...	13343116	2
Caelo	Prednisolone, microni...	13343105	2
Caelo	Prednisolone, microni...	13343106	2
Caelo	Prednisolone, microni...	15206102	1
Caelo	Prednisolone, microni...	13343105	2
Caelo	Prednisolone, microni...	121202102	2
Caelo	Prednisolone, microni...	13343106	5
Caelo	Prednisolone, microni...	15400003	21
Caelo	Prednisolone, microni...	9193501	1
Caelo	Prednisolone, microni...	15400002	17
Caelo	Prednisolone, microni...	15206106	1
Caelo	Prednisolone, microni...	k04c20121219	1
Caelo	Prednisolone, microni...	15206108	12
Caelo	Prednisolone, microni...	13343114	5
Caelo	Prednisolone, microni...	133414	1
Caelo	Prednisolone, microni...	1410B03302265	1
Caelo	Prednisolone, microni...	13343101	1
Caelo	Prednisolone, microni...	10196606	1
Caelo	Prednisolone, microni...	15206101	9
Caelo	Prednisolone, microni...	13343104	8
Caelo	Prednisolone, microni...	13343113	1
Caelo	Prednisolone, microni...	13343106	1
Caelo	Prednisolone, microni...	15206101	1
Caelo	Prednisolone, microni...	10196626	1
Caelo	Prednisolone, microni...	15206102	16
Caelo	Prednisolone, microni...	11/26-N04	1
Caelo	Prednisolone, microni...	15206105	7
Caelo	Prednisolone, microni...	15206102	1
Caelo	Prednisolone, microni...	13343111	5
Caelo	Prednisolone, microni...	13343115	7
Caelo	Prednisolone, microni...	16034301	2
Caelo	Prednisolone, microni...	10196627	7
Caelo	Prednisolone, microni...	12324221	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Prednisolone, microni...	2413Q02105	1
Caelo	Prednisolone, microni...	12324207	4
Caelo	Prednisolone, microni...	1234205	1
Caelo	Prednisolone, microni...	16134315, 180508P	1
Caelo	Prednisolone, microni...	17f07-b05	1
Caelo	Prednisolone, microni...	16134305	12
Caelo	Prednisolone, microni...	15400003	1
Caelo	Prednisolone, microni...	16134315	8
Caelo	Prednisolone, microni...	12324213	8
Caelo	Prednisolone, microni...	702020-01	1
Caelo	Prednisolone, microni...	17277703	5
Caelo	Prednisolone, microni...	17277704	9
Caelo	Prednisolone, microni...	NEPL180603	1
Caelo	Prednisolone, microni...	15055101	6
Caelo	Prednisolone, microni...	17277706	1
Caelo	Prednisolone, microni...	15055102	1
Caelo	Prednisolone, microni...	15055102	4
Caelo	Prednisolone, microni...	10196611	2
Caelo	Prednisolone, microni...	12324205	13
Caelo	Prednisolone, microni...	17277702	3
Caelo	Prednisolone, microni...	17020-01	1
Caelo	Prednisolone, microni...	12324226	3
Caelo	Prednisolone, microni...	12324220	1
Caelo	Prednisolone, microni...	15400010	4
Caelo	Prednisolone, microni...	12324221	2
Caelo	Prednisolone, microni...	1232417	2
Caelo	Prednisolone, microni...	12324219	9
Caelo	Prednisolone, microni...	16134306	4
Caelo	Prednisolone, microni...	161343032	1
Caelo	Prednisolone, microni...	411092-1-P	1
Caelo	Prednisolone, microni...	16134310	13
Caelo	Prednisolone, microni...	16134303	1
Caelo	Prednisolone, microni...	14G30-B14-302262	5
Caelo	Prednisolone, microni...	18138113	1
Caelo	Prednisolone, microni...	1613414	1
Caelo	Prednisolone, microni...	16134314	9
Caelo	Prednisolone, microni...	16134316	12
Caelo	Prednisolone, microni...	12324203	3
Caelo	Prednisolone, microni...	1511011	1
Caelo	Prednisolone, microni...	4513U2105	1
Caelo	Prednisolone, microni...	10196604	1
Caelo	Prednisolone, microni...	12324226	3
Caelo	Prednisolone, microni...		1
Caelo	Prednisolone, microni...	12324217	6
Caelo	Prednisolone, microni...	12324204	4
Caelo	Prednisolone, microni...	20021303	1
Caelo	Prednisolone, microni...	13343116	2
Caelo	Prednisolone, microni...	15206103	10
Caelo	Prednisolone, microni...	10196626	7
Caelo	Prednisolone, microni...	12324202	9
Caelo	Prednisolone, microni...	12185702	2
Caelo	Prednisolone, microni...	NEPL160202	1
Caelo	Prednisolone, microni...	13343115	1
Caelo	Prednisolone, microni...	60367086	2
Caelo	Prednisolone, microni...	12171301	6
Caelo	Prednisolone, microni...	12171306	2
Caelo	Prednisolone, microni...	12171302	1
Caelo	Prednisolone, microni...	15400001	2
Caelo	Prednisolone, microni...	16134302	12
Caelo	Prednisolone, microni...	4712I-02105	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Prednisolone, microni...	16C15-B01_320761	1
Caelo	Prednisolone, microni...	17277716	1
Caelo	Prednisolone, microni...	12171303	2
Caelo	Prednisolone, microni...	16C15-B01-320761	5
Caelo	Prednisolone, microni...	1727710	1
Caelo	Prednisolone, microni...	1540003	1
Caelo	Prednisolone, microni...	10196625	1
Caelo	Prednisolone, microni...	10196622	2
Caelo	Prednisolone, microni...	17277711	2
Caelo	Prednisolone, microni...	15400004	16
Caelo	Prednisone	204312107	1
Caelo	Prednisone	12177604	1
Caelo	Prednisone	16277206	2
Caelo	Prednisone	17166504	1
Caelo	Prednisone	17166507	3
Caelo	Prednisone	10037817	1
Caelo	Prednisone	16277201	1
Caelo	Prednisone	15218607	1
Caelo	Prednisone	17166502	1
Caelo	Prednisone	15218606	3
Caelo	Prednisone	5087	1
Caelo	Prednisone	13262304	3
Caelo	Prednisone	2114Q2107	1
Caelo	Prednisone	12232705	3
Caelo	Prednisone	12232713	1
Caelo	Prednisone	15218602	4
Caelo	Prednisone	10037826	1
Caelo	Prilocaine hydrochlor...	15230304	1
Caelo	Prilocaine hydrochlor...	14A03-B34-290281	1
Caelo	Primrose flower	11374404	1
Caelo	Primrose flower	13086603	1
Caelo	Primrose flower	11206808	1
Caelo	Primrose flower	5021303	1
Caelo	Primrose flower	11374403	1
Caelo	Primrose flower	13086608	1
Caelo	Primrose flower	13276904	1
Caelo	Primrose flower	9021301	1
Caelo	Primrose flower	1308660	1
Caelo	Primrose flower	11374402	1
Caelo	Primrose flower	25091306	1
Caelo	Primrose flower	10207104	1
Caelo	Primrose flower	10072509	1
Caelo	Primrose flower	11206807	1
Caelo	Primrose flower	13086603	5
Caelo	Primrose flower	13086605	1
Caelo	Primrose flower	11206804	1
Caelo	Primrose flower	11206811	1
Caelo	Primrose flower	113744044	1
Caelo	Primrose flower	11206810	5
Caelo	Primrose flower	11206803	4
Caelo	Primrose flower	11374401	10
Caelo	Primrose flower	13086602	5
Caelo	Primrose flower	11374403	3
Caelo	Primrose flower	11206811	4
Caelo	Primrose flower	13086602	3
Caelo	Primrose root, cut	12064506	4
Caelo	Primrose root, cut	12064503	1
Caelo	Primrose root, cut	13087001	1
Caelo	Primrose root, cut	4021410	1
Caelo	Primrose root, cut	10382010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Primrose root, cut	10382009	1
Caelo	Primrose root, cut	10382013	1
Caelo	Primrose root, cut	24011301	1
Caelo	Primrose root, cut	13087005-3534411072	1
Caelo	Primrose root, cut	10382012	1
Caelo	Primrose root, cut	Int.Ch:318090	1
Caelo	Primrose root, cut	13087004	1
Caelo	Primrose root, cut	13067001	1
Caelo	Primrose root, cut	71653517	1
Caelo	Primrose root, cut	5010E-02159	1
Caelo	Primrose root, cut	13087005	3
Caelo	Primrose root, cut	13087002	4
Caelo	Procaine hydrochloride	17343311	1
Caelo	Procaine hydrochloride	17167505	1
Caelo	Procaine hydrochloride	18041806	1
Caelo	Procaine hydrochloride	17167504	1
Caelo	Procaine hydrochloride	14164409	1
Caelo	Procaine hydrochloride	16137812	3
Caelo	Procaine hydrochloride	13130-N06	1
Caelo	Procaine hydrochloride	15006912	3
Caelo	Procaine hydrochloride	15006913	1
Caelo	Procaine hydrochloride	15006909	1
Caelo	Procaine hydrochloride	13F12-B40-302510	1
Caelo	Procaine hydrochloride	15006906	4
Caelo	Procaine hydrochloride	10159723	1
Caelo	Procaine hydrochloride	16173906	1
Caelo	Procaine hydrochloride	10159729	1
Caelo	Procaine hydrochloride	14164408	1
Caelo	Procaine hydrochloride	11370204	1
Caelo	Procaine hydrochloride	12177801	6
Caelo	Procaine hydrochloride	15006902	1
Caelo	Procaine hydrochloride	11370206	1
Caelo	Procaine hydrochloride	10159721	1
Caelo	Procaine hydrochloride	11370215	1
Caelo	Procaine hydrochloride	11370219	1
Caelo	Procaine hydrochloride	131099903	1
Caelo	Procaine hydrochloride	11370216	2
Caelo	Procaine hydrochloride	13109908	2
Caelo	Procaine hydrochloride	13109915	1
Caelo	Procaine hydrochloride	13109906	1
Caelo	Procaine hydrochloride	13/09-N08	1
Caelo	Procaine hydrochloride	13109905	1
Caelo	Procaine hydrochloride	7011506	1
Caelo	Procaine hydrochloride	13109916	2
Caelo	Procaine hydrochloride	60998127	2
Caelo	Procaine hydrochloride	11370203	1
Caelo	Progesterone, microni...	18098801	25
Caelo	Progesterone, microni...	2006500002	1
Caelo	Progesterone, microni...	15I11-B05-314184	6
Caelo	Progesterone, microni...	1511-B05-313598	1
Caelo	Progesterone, microni...	15319601	5
Caelo	Progesterone, microni...	15124910	3
Caelo	Progesterone, microni...	13HI01-N03	1
Caelo	Progesterone, microni...	14253929	1
Caelo	Progesterone, microni...	12195924	3
Caelo	Progesterone, microni...	14253915	1
Caelo	Progesterone, microni...	14253911	2
Caelo	Progesterone, microni...	17348407	2
Caelo	Progesterone, microni...	16103604	1
Caelo	Progesterone, microni...	13155929	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Progesterone, microni...	13155928	3
Caelo	Progesterone, microni...	13155920	2
Caelo	Progesterone, microni...	14H13-B05-3002016	1
Caelo	Progesterone, microni...	15C31-808	1
Caelo	Progesterone, microni...	160100902	1
Caelo	Progesterone, microni...	15319602	1
Caelo	Progesterone, microni...	16103606	9
Caelo	Progesterone, microni...	16103602	2
Caelo	Progesterone, microni...	16015B14-315279	1
Caelo	Progesterone, microni...	16103610	4
Caelo	Progesterone, microni...	16103613	20
Caelo	Progesterone, microni...	13155927	2
Caelo	Progesterone, microni...	15319618	2
Caelo	Progesterone, microni...	15319608	6
Caelo	Progesterone, microni...	16D13-B04	1
Caelo	Progesterone, microni...	16103601	2
Caelo	Progesterone, microni...	15319607	1
Caelo	Progesterone, microni...	16103609	30
Caelo	Progesterone, microni...	161036303	9
Caelo	Progesterone, microni...	16A20-B02-319103	3
Caelo	Progesterone, microni...	15124919	2
Caelo	Progesterone, microni...	16103603	11
Caelo	Progesterone, microni...	258/61	1
Caelo	Progesterone, microni...	4814A-03235	1
Caelo	Progesterone, microni...	18f22-B04-354899	1
Caelo	Progesterone, microni...	17348402	1
Caelo	Progesterone, microni...	15319619	1
Caelo	Progesterone, microni...	12195913	3
Caelo	Progesterone, microni...	11177221	2
Caelo	Progesterone, microni...	17/21-B03-343284	2
Caelo	Progesterone, microni...	16L29-B03-330995	2
Caelo	Progesterone, microni...	16103607	2
Caelo	Progesterone, microni...	EU-HTT160708	1
Caelo	Progesterone, microni...	17149906	10
Caelo	Progesterone, microni...	17060907	10
Caelo	Progesterone, microni...	17149907	1
Caelo	Progesterone, microni...	17149904	9
Caelo	Progesterone, microni...	B03333671	1
Caelo	Progesterone, microni...	16J05-B04-328263	5
Caelo	Progesterone, microni...	1719907	10
Caelo	Progesterone, microni...	17060903	5
Caelo	Progesterone, microni...	4817Q-03235	1
Caelo	Progesterone, microni...	1760911	1
Caelo	Progesterone, microni...	1706093	1
Caelo	Progesterone, microni...	15319620	1
Caelo	Progesterone, microni...	17060909	6
Caelo	Progesterone, microni...	17277801	22
Caelo	Progesterone, microni...	17060905	1
Caelo	Progesterone, microni...	1721bo3342134	1
Caelo	Progesterone, microni...	17149903	1
Caelo	Progesterone, microni...	9195408	1
Caelo	Progesterone, microni...	15319603	42
Caelo	Progesterone, microni...	9195429	1
Caelo	Progesterone, microni...	13155921	3
Caelo	Progesterone, microni...	12195828	1
Caelo	Progesterone, microni...	140100201	1
Caelo	Progesterone, microni...	17060901	41
Caelo	Progesterone, microni...	13155902	11
Caelo	Progesterone, microni...	12195912	4
Caelo	Progesterone, microni...	130404002	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Progesterone, microni...	17348404	5
Caelo	Progesterone, microni...	C01-120613	1
Caelo	Progesterone, microni...	12195928	8
Caelo	Progesterone, microni...	12195914	3
Caelo	Progesterone, microni...	12195920	8
Caelo	Progesterone, microni...	17060904	1
Caelo	Progesterone, microni...	12195906	4
Caelo	Progesterone, microni...	11177219	3
Caelo	Progesterone, microni...	12195905	3
Caelo	Progesterone, microni...	18B06-B10	3
Caelo	Progesterone, microni...	11177207	1
Caelo	Progesterone, microni...	12195926	19
Caelo	Promethazine hydrochl...	11352305	1
Caelo	Promethazine hydrochl...	13119211	1
Caelo	Promethazine hydrochl...	17058601	1
Caelo	Promethazine hydrochl...	19051805	1
Caelo	Propranolol hydrochlo...	16138302	7
Caelo	Propranolol hydrochlo...	15283302	7
Caelo	Propranolol hydrochlo...	15F16-B02	1
Caelo	Propranolol hydrochlo...	173330015654	1
Caelo	Propranolol hydrochlo...	1542251	1
Caelo	Propranolol hydrochlo...	15422501	12
Caelo	Propranolol hydrochlo...	15422502	15
Caelo	Propranolol hydrochlo...	154225	1
Caelo	Propranolol hydrochlo...	17K22-B03-345013	3
Caelo	Propranolol hydrochlo...	15422505	7
Caelo	Propranolol hydrochlo...	1733301	1
Caelo	Propranolol hydrochlo...	M141203	2
Caelo	Propranolol hydrochlo...	17K22-B03-344247	2
Caelo	Propranolol hydrochlo...	17333005	1
Caelo	Propranolol hydrochlo...	15k11-B01	1
Caelo	Propranolol hydrochlo...	17333001	13
Caelo	Propranolol hydrochlo...	16E09-B09-330547	8
Caelo	Propranolol hydrochlo...	150082405	1
Caelo	Propranolol hydrochlo...	17172203	12
Caelo	Propranolol hydrochlo...	16138306	2
Caelo	Propranolol hydrochlo...	17172204	11
Caelo	Propranolol hydrochlo...	15F16-B02-312979	7
Caelo	Propranolol hydrochlo...	17172205	3
Caelo	Propranolol hydrochlo...	15283301	11
Caelo	Propranolol hydrochlo...	15082404	1
Caelo	Propranolol hydrochlo...	18179403	1
Caelo	Propranolol hydrochlo...	18179402	16
Caelo	Propranolol hydrochlo...	17K22-B03-347251	3
Caelo	Propranolol hydrochlo...	17333006	4
Caelo	Propranolol hydrochlo...	12a27-n04	1
Caelo	Propranolol hydrochlo...	12370703	1
Caelo	Propranolol hydrochlo...	M130104	1
Caelo	Propranolol hydrochlo...	15082405	2
Caelo	Propranolol hydrochlo...	12370709	1
Caelo	Propranolol hydrochlo...	12370710	1
Caelo	Propranolol hydrochlo...	m130104	1
Caelo	Propranolol hydrochlo...	12370712	1
Caelo	Propranolol hydrochlo...	13379804	1
Caelo	Propranolol hydrochlo...	13/13-N04	1
Caelo	Propranolol hydrochlo...	13379806	1
Caelo	Propranolol hydrochlo...	13379807	1
Caelo	Propranolol hydrochlo...	13379809	2
Caelo	Propranolol hydrochlo...	13379809	1
Caelo	Propranolol hydrochlo...	13I13N04	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Propranolol hydrochlo...	1712204	1
Caelo	Propyl-4-hydroxybenzo...	13051504	1
Caelo	Propyl-4-hydroxybenzo...	11373524	1
Caelo	Propyl-4-hydroxybenzo...	11373518	2
Caelo	Propyl-4-hydroxybenzo...	82792159	1
Caelo	Propyl-4-hydroxybenzo...	13365009	3
Caelo	Propyl-4-hydroxybenzo...	13365006	2
Caelo	Propyl-4-hydroxybenzo...	13365006	3
Caelo	Propyl-4-hydroxybenzo...	113735223	1
Caelo	Propyl-4-hydroxybenzo...	221014PAS01	1
Caelo	Propyl-4-hydroxybenzo...	3041408	1
Caelo	Propyl-4-hydroxybenzo...	13100217	2
Caelo	Propyl-4-hydroxybenzo...	13100218	2
Caelo	Propyl-4-hydroxybenzo...	13100208	1
Caelo	Propyl-4-hydroxybenzo...	11102609	1
Caelo	Propyl-4-hydroxybenzo...	1913I-01931	1
Caelo	Propyl-4-hydroxybenzo...	13365011	4
Caelo	Propyl-4-hydroxybenzo...	713Q-01931	1
Caelo	Propyl-4-hydroxybenzo...	290114PML03	1
Caelo	Propyl-4-hydroxybenzo...	13365011	1
Caelo	Propyl-4-hydroxybenzo...	11373505	2
Caelo	Propyl-4-hydroxybenzo...	15246502	2
Caelo	Propyl-4-hydroxybenzo...	14359019	1
Caelo	Propyl-4-hydroxybenzo...	3115M-01931	1
Caelo	Propyl-4-hydroxybenzo...	15246501	1
Caelo	Propyl-4-hydroxybenzo...	14359006	3
Caelo	Propyl-4-hydroxybenzo...	16359006	1
Caelo	Propyl-4-hydroxybenzo...	14359005	9
Caelo	Propyl-4-hydroxybenzo...	13365012	1
Caelo	Propyl-4-hydroxybenzo...	14359009	1
Caelo	Propyl-4-hydroxybenzo...	4314U-01931	1
Caelo	Propyl-4-hydroxybenzo...	13100215	1
Caelo	Propyl-4-hydroxybenzo...	11102611	2
Caelo	Propyl-4-hydroxybenzo...	13100219	1
Caelo	Propyl-4-hydroxybenzo...	14359008	22
Caelo	Propyl-4-hydroxybenzo...	14359020	5
Caelo	Propyl-4-hydroxybenzo...	15276514	2
Caelo	Propyl-4-hydroxybenzo...	16259505	7
Caelo	Propyl-4-hydroxybenzo...	13100215	9
Caelo	Propyl-4-hydroxybenzo...	11102624	1
Caelo	Propyl-4-hydroxybenzo...	13365012	1
Caelo	Propyl-4-hydroxybenzo...	17343108	1
Caelo	Propyl-4-hydroxybenzo...	17343105	1
Caelo	Propyl-4-hydroxybenzo...	17343110	1
Caelo	Propyl-4-hydroxybenzo...	17343101	5
Caelo	Propyl-4-hydroxybenzo...	113735	1
Caelo	Propyl-4-hydroxybenzo...	11121709	1
Caelo	Propyl-4-hydroxybenzo...	16259516	12
Caelo	Propyl-4-hydroxybenzo...	1108211931	1
Caelo	Propyphenazone	12039101	2
Caelo	Propyphenazone	12039103	2
Caelo	Propyphenazone	12039112	1
Caelo	Propyphenazone	13036702	2
Caelo	Propyphenazone	13332205	3
Caelo	Propyphenazone	13137104	1
Caelo	Propyphenazone	15122502	1
Caelo	Propyphenazone	17351305	1
Caelo	Propyphenazone	15348913	1
Caelo	Propyphenazone	15348912	2
Caelo	Propyphenazone	15348908	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Propyphenazone	15348904	2
Caelo	Propyphenazone	15348902	1
Caelo	Propyphenazone	15122504	2
Caelo	Propyphenazone	15122505	1
Caelo	Psyllium seed	13305903	1
Caelo	Psyllium seed	11185803	1
Caelo	Psyllium seed	11185811	4
Caelo	Psyllium seed	13305904	2
Caelo	Psyllium seed	10291205	1
Caelo	Psyllium seed	12140301	7
Caelo	Psyllium seed	10291206	1
Caelo	Psyllium seed	12185504	10
Caelo	Psyllium seed	10291207	6
Caelo	Psyllium seed	10291204	1
Caelo	Psyllium seed	13234309	1
Caelo	Psyllium seed	13305910	4
Caelo	Psyllium seed	12185503	3
Caelo	Psyllium seed	13305907	5
Caelo	Psyllium seed	15190303	4
Caelo	Psyllium seed	15190304	3
Caelo	Psyllium seed	10061503	1
Caelo	Psyllium seed	7380a120911	1
Caelo	Psyllium seed	13234310	3
Caelo	Psyllium seed	12276510	1
Caelo	Psyllium seed	12185502	1
Caelo	Psyllium seed	12185508	4
Caelo	Psyllium seed	12185509	5
Caelo	Psyllium seed	15190306	1
Caelo	Psyllium seed	13305906	2
Caelo	Pumpkin seed, peeled	16081206	2
Caelo	Pumpkin seed, peeled	16081202	1
Caelo	Pumpkin seed, peeled	14186910	1
Caelo	Pumpkin seed, peeled	12249001	2
Caelo	Pumpkin seed, peeled	12249011	2
Caelo	Pumpkin seed, peeled	122490911	1
Caelo	Pumpkin seed, peeled	140736306	1
Caelo	Pumpkin seed, peeled	14186905	1
Caelo	Pumpkin seed, peeled	14073603	1
Caelo	Pyridoxine hydrochlor...	13102510	2
Caelo	Pyridoxine hydrochlor...	16222506	2
Caelo	Pyridoxine hydrochlor...	16222521	1
Caelo	Pyridoxine hydrochlor...	16222520	1
Caelo	Pyridoxine hydrochlor...	16222514	2
Caelo	Pyridoxine hydrochlor...	16222501	5
Caelo	Pyridoxine hydrochlor...	15218709	1
Caelo	Pyridoxine hydrochlor...	11021406	1
Caelo	Pyridoxine hydrochlor...	15218713	3
Caelo	Pyridoxine hydrochlor...	15218710	2
Caelo	Pyridoxine hydrochlor...	15218702	6
Caelo	Pyridoxine hydrochlor...	14262306	1
Caelo	Pyridoxine hydrochlor...	13102516	1
Caelo	Pyridoxine hydrochlor...	12102510	1
Caelo	Pyridoxine hydrochlor...	13102509	2
Caelo	Pyridoxine hydrochlor...	15218707	1
Caelo	Pyridoxine hydrochlor...	16222515	2
Caelo	Pyridoxine hydrochlor...	14262302	3
Caelo	Quinine hydrochloride	16283305	1
Caelo	Quinine hydrochloride	153776201	1
Caelo	Quinine hydrochloride	12316611	1
Caelo	Quinine hydrochloride	13D24-N06	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Quinine hydrochloride	12316605	1
Caelo	Quinine hydrochloride	12316604	1
Caelo	Quinine hydrochloride	16283302	1
Caelo	Quinine sulfate dihyd...	15B06B03304359	1
Caelo	Quinoline yellow	14274305	1
Caelo	Quinoline yellow	17034104	6
Caelo	Quinoline yellow	11089606	1
Caelo	Quinoline yellow	17034102	5
Caelo	Quinoline yellow	16050106	3
Caelo	Quinoline yellow	713E-03146	1
Caelo	Quinoline yellow	15276801	2
Caelo	Quinoline yellow	12042007	3
Caelo	Quinoline yellow	20301211-6	1
Caelo	Quinoline yellow	14274309	1
Caelo	Quinoline yellow	12042010	1
Caelo	Quinoline yellow	15276804	2
Caelo	Quinoline yellow	15276802	1
Caelo	Quinoline yellow	14363202	3
Caelo	Quinoline yellow	16050101	5
Caelo	Quinoline yellow	11089610	1
Caelo	Quinoline yellow	118E-03146	1
Caelo	Quinoline yellow	14363207	1
Caelo	Quinoline yellow	17227706	2
Caelo	Quinoline yellow	13109102	2
Caelo	Quinoline yellow	13109109	1
Caelo	Quinoline yellow	17227702	1
Caelo	Quinoline yellow	12042002	3
Caelo	Quinoline yellow	16050108	1
Caelo	Quinoline yellow	14049401	2
Caelo	Quinoline yellow	17227703	3
Caelo	Quinoline yellow	17227701	5
Caelo	Quinoline yellow	17034101	1
Caelo	Quinoline yellow	11089608	1
Caelo	Quinolinol sulfate po...	15147714	1
Caelo	Quinolinol sulfate po...	16221811	1
Caelo	Quinolinol sulfate po...	16221809	3
Caelo	Quinolinol sulfate po...	16221807	2
Caelo	Quinolinol sulfate po...	15147708	4
Caelo	Quinolinol sulfate po...	15147706	5
Caelo	Quinolinol sulfate po...	9106304	1
Caelo	Quinolinol sulfate po...	10281601	1
Caelo	Quinolinol sulfate po...	15147704	4
Caelo	Quinolinol sulfate po...	13403501	3
Caelo	Quinolinol sulfate po...	7051305	1
Caelo	Quinolinol sulfate po...	13071101	2
Caelo	Quinolinol sulfate po...	13071104	4
Caelo	Quinolinol sulfate po...	13403513	1
Caelo	Quinolinol sulfate po...	13403501	4
Caelo	Quinolinol sulfate po...	13071106	1
Caelo	Quinolinol sulfate po...	15147701	3
Caelo	Quinolinol sulfate po...	13403506	12
Caelo	Ramsons herb, cut	11182204	1
Caelo	Ramsons herb, cut	12380803	1
Caelo	Ramsons herb, cut	11182207	1
Caelo	Ramsons herb, cut	12380801	2
Caelo	Ramsons herb, cut	12380802	2
Caelo	Ramsons herb, cut	11182203	1
Caelo	Ramsons herb, cut	14343305	1
Caelo	Raspberry leaf, cut	12103111	15
Caelo	Raspberry leaf, cut	12130107	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Raspberry leaf, cut	12362702	7
Caelo	Raspberry leaf, cut	8031307	1
Caelo	Raspberry leaf, cut	18041304	1
Caelo	Raspberry leaf, cut	25091201	1
Caelo	Raspberry leaf, cut	1210315	1
Caelo	Raspberry leaf, cut	713E-01511	1
Caelo	Raspberry leaf, cut	2909211511	1
Caelo	Raspberry leaf, cut	12362706	28
Caelo	Raspberry leaf, cut	10013410	1
Caelo	Raspberry leaf, cut	17121203	1
Caelo	Raspberry leaf, cut	11121201	1
Caelo	Raspberry leaf, cut	22081201	1
Caelo	Raspberry leaf, cut	9101204	1
Caelo	Raspberry leaf, cut	9081205	1
Caelo	Raspberry leaf, cut	12103101	7
Caelo	Raspberry leaf, cut	18101203	1
Caelo	Raspberry leaf, cut	12103108	8
Caelo	Raspberry leaf, cut	10013411	1
Caelo	Raspberry leaf, cut	307211511	1
Caelo	Raspberry leaf, cut	10013409	2
Caelo	Raspberry leaf, cut	10013407	1
Caelo	Raspberry leaf, cut	11030706	3
Caelo	Raspberry leaf, cut	12103109	2
Caelo	Raspberry leaf, cut	21051302	1
Caelo	Raspberry leaf, cut	2271C120508	5
Caelo	Raspberry leaf, cut	13127901	3
Caelo	Raspberry leaf, cut	314U1511	1
Caelo	Raspberry leaf, cut	1310603	1
Caelo	Raspberry leaf, cut	9011404	1
Caelo	Raspberry leaf, cut	13038002	1
Caelo	Raspberry leaf, cut	7011402	1
Caelo	Raspberry leaf, cut	13101602-347312141	1
Caelo	Raspberry leaf, cut	13101604	3
Caelo	Raspberry leaf, cut	12362708	4
Caelo	Raspberry leaf, cut	12362711	9
Caelo	Raspberry leaf, cut	12362710-3432311160	1
Caelo	Raspberry leaf, cut	12362709	7
Caelo	Raspberry leaf, cut	231013c	1
Caelo	Raspberry leaf, cut	13101803	1
Caelo	Raspberry leaf, cut	13038003	1
Caelo	Raspberry leaf, cut	12362710	13
Caelo	Raspberry leaf, cut	10136404	1
Caelo	Raspberry leaf, cut	12392709	1
Caelo	Raspberry leaf, cut	10136412	1
Caelo	Raspberry leaf, cut	11071301	1
Caelo	Raspberry leaf, cut	11030704	2
Caelo	Raspberry leaf, cut	IntCh:311090	1
Caelo	Raspberry leaf, cut	1236705	1
Caelo	Raspberry leaf, cut	12362707	4
Caelo	Raspberry leaf, cut	11030705	1
Caelo	Raspberry leaf, cut	12362704	5
Caelo	Raspberry leaf, cut	12368706	1
Caelo	Raspberry leaf, cut	13101603	5
Caelo	Raspberry leaf, cut	20081311	1
Caelo	Raspberry leaf, cut	14018803	1
Caelo	Raspberry leaf, cut	6859427	1
Caelo	Raspberry leaf, cut	13101605	1
Caelo	Raspberry leaf, cut	30011408	1
Caelo	Raspberry leaf, cut	12362705	9
Caelo	Raspberry leaf, cut	11030707	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Red poppy petals, cut	12330205	1
Caelo	Red-rose petal	12380004	1
Caelo	Red-rose petal	12380008	2
Caelo	Red-rose petal, cut	13194402	1
Caelo	Red-rose petal, cut	13194403	4
Caelo	Red-rose petal, cut	12380008	1
Caelo	Red-rose petal, cut	10113211	1
Caelo	Red-rose petal, cut	10306514	1
Caelo	Red-rose petal, cut	11304901	3
Caelo	Red-rose petal, cut	12380007-3357319071	1
Caelo	Red-rose petal, cut	11304904	4
Caelo	Red-rose petal, cut	11011301	1
Caelo	Red-rose petal, cut	12380001	4
Caelo	Red-rose petal, cut	11304906	1
Caelo	Red-rose petal, cut	12380004	2
Caelo	Red-rose petal, cut	11304912	2
Caelo	Redwood	12041403	1
Caelo	Redwood	13106901	1
Caelo	Redwood	12049607	2
Caelo	Redwood	12049608	2
Caelo	Resocrinol	10218809	2
Caelo	Resocrinol	13293110	1
Caelo	Resocrinol	12255614	2
Caelo	Resocrinol	19061305	1
Caelo	Resocrinol	206370-0003	1
Caelo	Resocrinol	13293115	2
Caelo	Resocrinol	15009502	1
Caelo	Resocrinol	15009507	6
Caelo	Resocrinol	15009511	4
Caelo	Resocrinol	18111703	1
Caelo	Resocrinol	13293103	1
Caelo	Resocrinol	16D29-B01	2
Caelo	Resocrinol	16153307	1
Caelo	Resocrinol	18015601	1
Caelo	Resocrinol	16153312	1
Caelo	Resocrinol	11309610	4
Caelo	Resocrinol	18216601	1
Caelo	Resocrinol	16153309	6
Caelo	Resocrinol	16271807	1
Caelo	Resocrinol	18130901	1
Caelo	Resocrinol	10111210	1
Caelo	Resocrinol	2041512	1
Caelo	Resocrinol	11d29-n14	2
Caelo	Resocrinol	5052	1
Caelo	Resocrinol	15009515	1
Caelo	Resocrinol	13/18-B19-291534	1
Caelo	Resocrinol	18216609	1
Caelo	Resocrinol	18216603	1
Caelo	Restharrow root, cut	11247302	2
Caelo	Restharrow root, cut	15101313	1
Caelo	Restharrow root, cut	12330003	2
Caelo	Restharrow root, cut	12330005	5
Caelo	Restharrow root, cut	72439069	1
Caelo	Restharrow root, cut	71035467	1
Caelo	Restharrow root, cut	Int.Ch:212150	1
Caelo	Restharrow root, cut	6101201	1
Caelo	Restharrow root, cut	12330002	1
Caelo	Restharrow root, cut	11103107	2
Caelo	Restharrow root, cut	20111213	1
Caelo	Restharrow root, cut	21011306	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Restharrow root, cut	12330001	14
Caelo	Restharrow root, cut	1233001	2
Caelo	Restharrow root, cut	11103105	2
Caelo	Restharrow root, cut	4021406	1
Caelo	Restharrow root, cut	1230004	1
Caelo	Rhubarb, cut	11057804	1
Caelo	Rhubarb, cut	10071302	1
Caelo	Rhubarb, cut	12211902	3
Caelo	Rhubarb, cut	11057810	1
Caelo	Rhubarb, cut	11321806	1
Caelo	Rhubarb, cut	25111305	1
Caelo	Rhubarb, cut	12211903	4
Caelo	Rhubarb, cut	13235201	2
Caelo	Rhubarb, cut	13235203	1
Caelo	Rhubarb, cut	11057818	2
Caelo	Rhubarb, cut	12211904	2
Caelo	Rhubarb, cut	13235205	1
Caelo	Rhubarb, cut	12211906	1
Caelo	Rhubarb, cut	15101301	1
Caelo	Rhubarb, cut	11265902	1
Caelo	Rhubarb, cut	14155404	1
Caelo	Rhubarb, cut	12211908	1
Caelo	Rhubarb, powdered	13193106	2
Caelo	Rhubarb, powdered	23061402	1
Caelo	Rhubarb, powdered	13193109	1
Caelo	Rhubarb, powdered	13193104	4
Caelo	Rhubarb, powdered	11043825	2
Caelo	Rhubarb, powdered	11043819	2
Caelo	Rhubarb, powdered	11043816	1
Caelo	Rhubarb, powdered	11043822	1
Caelo	Rhubarb, powdered	9071201	1
Caelo	Rhubarb, powdered	11043816	2
Caelo	Rhubarb, powdered	11043817	1
Caelo	Rhubarb, powdered	11043818	1
Caelo	Rhubarb, powdered	11043818	1
Caelo	Rhubarb, powdered	11043826	1
Caelo	Riboflavin	17G14-B05-337885	1
Caelo	Riboflavin	17272204	1
Caelo	Riboflavin	15330708	1
Caelo	Ribwort plantain, cut	13113006	6
Caelo	Ribwort plantain, cut	12362906	7
Caelo	Ribwort plantain, cut	4021411	1
Caelo	Ribwort plantain, cut	1313007	1
Caelo	Ribwort plantain, cut	16101304	1
Caelo	Ribwort plantain, cut	25041301	1
Caelo	Ribwort plantain, cut	1362903	1
Caelo	Ribwort plantain, cut	13113001	8
Caelo	Ribwort plantain, cut	11061304	1
Caelo	Ribwort plantain, cut	13113007	6
Caelo	Ribwort plantain, cut	15111302	1
Caelo	Ribwort plantain, cut	20111304	1
Caelo	Ribwort plantain, cut	13232401	1
Caelo	Ribwort plantain, cut	13113010	4
Caelo	Ribwort plantain, cut	28061301	1
Caelo	Ribwort plantain, cut	10011402	1
Caelo	Ribwort plantain, cut	13113004-3488312172	1
Caelo	Ribwort plantain, cut	23121304	1
Caelo	Ribwort plantain, cut	13113008	5
Caelo	Ribwort plantain, cut	25091305	1
Caelo	Ribwort plantain, cut	13113005	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Ribwort plantain, cut	12362901	1
Caelo	Ribwort plantain, cut	20111208	1
Caelo	Ribwort plantain, cut	12362903-3388311070	1
Caelo	Ribwort plantain, cut	10315715	8
Caelo	Ribwort plantain, cut	17071201	1
Caelo	Ribwort plantain, cut	10315711	4
Caelo	Ribwort plantain, cut	8061208	1
Caelo	Ribwort plantain, cut	10315714	2
Caelo	Ribwort plantain, cut	10315707	2
Caelo	Ribwort plantain, cut	10328702	2
Caelo	Ribwort plantain, cut	12225802	5
Caelo	Ribwort plantain, cut	15101204	1
Caelo	Ribwort plantain, cut	12225804	1
Caelo	Ribwort plantain, cut	19121202	1
Caelo	Ribwort plantain, cut	15096308	1
Caelo	Ribwort plantain, cut	15096301	1
Caelo	Ribwort plantain, cut	14199202	1
Caelo	Ribwort plantain, cut	10315705	1
Caelo	Ribwort plantain, cut	13410612	1
Caelo	Ribwort plantain, cut	12362902	7
Caelo	Ribwort plantain, cut	22021302	1
Caelo	Ribwort plantain, cut	13410601	1
Caelo	Ribwort plantain, cut	13410603	3
Caelo	Ribwort plantain, cut	2031302	1
Caelo	Ribwort plantain, cut	12225803	2
Caelo	Ribwort plantain, cut	10315712	2
Caelo	Rice starch	200613K3	1
Caelo	Rice starch	13154914	1
Caelo	Rice starch	13005905	2
Caelo	Rice starch	11132406	1
Caelo	Rice starch	10177507	1
Caelo	Rice starch	12214303	1
Caelo	Rice starch	12166007	3
Caelo	Rice starch	13154904	2
Caelo	Rice starch	9352108	1
Caelo	Rice starch	14025710	1
Caelo	Rice starch	12238714	1
Caelo	Rice starch	10061134	2
Caelo	Rice starch	10177515	1
Caelo	Rice starch	11186026	1
Caelo	Rice starch	14011612	3
Caelo	Rice starch	15119011	1
Caelo	Rice starch	13305108	1
Caelo	Rice starch	13005908	1
Caelo	Rice starch	2212M-01115	1
Caelo	Rice starch	16218002	1
Caelo	Rice starch	16144506	1
Caelo	Rice starch	16101401	1
Caelo	Rice starch	15118607	1
Caelo	Rice starch	17173501	2
Caelo	Rice starch	16218009	1
Caelo	Rice starch	17327503	2
Caelo	Rice starch	17327502	1
Caelo	Rice starch	14011630	3
Caelo	Rice starch	13305103	4
Caelo	Rice starch	12238719	1
Caelo	Rice starch	14011621	2
Caelo	Rice starch	14301214	1
Caelo	Rice starch	14025733	1
Caelo	Robert herb, cut	12319302	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Robert herb, cut	11147812	3
Caelo	Robert herb, cut	11147813	2
Caelo	Robert herb, cut	8121203	1
Caelo	Robert herb, cut	12319305	1
Caelo	Robert herb, cut	11147804	1
Caelo	Robert herb, cut	11147807	3
Caelo	Robert herb, cut	12319304	2
Caelo	Rose hips peel, cut	9121301	1
Caelo	Rose hips peel, cut	13317805	1
Caelo	Rose hips peel, cut	14050704	1
Caelo	Rose hips peel, cut	14050705	1
Caelo	Rose hips peel, cut	12357504	2
Caelo	Rose hips peel, cut	10157309	1
Caelo	Rose hips peel, cut	11301403	2
Caelo	Rose hips peel, cut	13317801	2
Caelo	Rose hips peel, cut	25061203	1
Caelo	Rose hips peel, cut	5091301	1
Caelo	Rose hips peel, cut	2851A120904konv.	1
Caelo	Rose hips peel, cut	12357505	6
Caelo	Rose hips peel, cut	12357502	2
Caelo	Rose hips peel, cut	12357501	1
Caelo	Rose hips peel, cut	12288301	1
Caelo	Rose hips peel, cut	15061210	1
Caelo	Rose hips peel, cut	12064207	1
Caelo	Rose hips peel, cut	7031305	1
Caelo	Rose hips peel, cut	11301407	12
Caelo	Rose hips peel, cut	10101206	1
Caelo	Rose hips peel, cut	11301405	2
Caelo	Rose hips peel, cut	12288302	6
Caelo	Rose hips peel, cut	11301401	1
Caelo	Rose hips with peel, ...	12064213	1
Caelo	Rose hips with peel, ...	19111204	1
Caelo	Rose hips with peel, ...	12064213	1
Caelo	Rose hips with peel, ...	114Q-01542	1
Caelo	Rose hips with peel, ...	12029008	1
Caelo	Rose hips with peel, ...	13097405	1
Caelo	Rose hips with peel, ...	13048703	1
Caelo	Rose hips with peel, ...	13048703	1
Caelo	Rose hips with peel, ...	13048701	2
Caelo	Rose hips with peel, ...	12064207	2
Caelo	Rose hips with peel, ...	12064206	1
Caelo	Rose hips with peel, ...	13048704	1
Caelo	Rose hips with peel, ...	12064209	1
Caelo	Rose hips with peel, ...	307211542	1
Caelo	Rose hips with peel, ...	12029011	1
Caelo	Rose hips with peel, ...	16111201	1
Caelo	Rose hips with peel, ...	12064209	4
Caelo	Rose hips with peel, ...	13048704	3
Caelo	Rose hips with peel, ...	1512211542	1
Caelo	Rose hips with peel, ...	12064207	3
Caelo	Rose hips with peel, ...	12064212	2
Caelo	Rose hips with peel, ...	12064210	3
Caelo	Rose hips with peel, ...	2103311542	1
Caelo	Rose hips with peel, ...	12064201	1
Caelo	Rose hips with peel, ...	12064212	3
Caelo	Rose hips with peel, ...	12064210	1
Caelo	Rose hips with peel, ...	12064202	1
Caelo	Rose hips with peel, ...	12064202	1
Caelo	Rose mallow flower, cut	13040401	2
Caelo	Rose mallow flower, cut	13040402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Rose mallow flower, cut	10244507	1
Caelo	Roselle	11036511	1
Caelo	Roselle	12007306	1
Caelo	Roselle	12267501	3
Caelo	Roselle	12267504	3
Caelo	Roselle	11036512	1
Caelo	Roselle	7121301	1
Caelo	Roselle	12267503	2
Caelo	Roselle	12267504-3428311150	1
Caelo	Roselle	12267502	3
Caelo	Roselle	11036504	1
Caelo	Roselle	12267504	1
Caelo	Roselle	12267502	4
Caelo	Roselle	11036509	1
Caelo	Roselle	12267507	2
Caelo	Rosemary leaf, cut	17352502	1
Caelo	Rosemary leaf, cut	15125008	1
Caelo	Rosemary leaf, cut	12225902	1
Caelo	Rosemary leaf, cut	1415901	1
Caelo	Rosemary leaf, cut	13128801	1
Caelo	Rosemary leaf, cut	28051404	1
Caelo	Rosemary leaf, cut	11193402	2
Caelo	Rosemary leaf, cut	11193405	1
Caelo	Rosemary leaf, cut	11193405	1
Caelo	Rosemary leaf, cut	14082901	1
Caelo	Rosemary leaf, cut	13128805	1
Caelo	Rosemary leaf, cut	13128806	2
Caelo	Rosemary leaf, cut	13128802	2
Caelo	Rosemary leaf, cut	11193404	1
Caelo	Rosemary leaf, cut	15125001	2
Caelo	Rosemary leaf, cut	14159601	1
Caelo	Rosemary leaf, cut	5041302	1
Caelo	Rosemary leaf, cut	11193406	1
Caelo	Rosemary leaf, cut	IntCh:313060	1
Caelo	Rosemary leaf, cut	11193408	6
Caelo	Rue herb, cut	11371803	1
Caelo	Rue herb, cut	12148801	1
Caelo	Rue herb, cut	12148806	1
Caelo	Rue herb, cut	12148804	1
Caelo	Rue herb, cut	12148803	2
Caelo	Rue herb, cut	9336802	1
Caelo	Rue herb, cut		1
Caelo	Saccharin sodium	11172805	2
Caelo	Saccharin sodium	90542309	1
Caelo	Saccharin sodium	11172814	1
Caelo	Saccharin sodium	1236207	1
Caelo	Saccharin sodium	11025605	1
Caelo	Saccharin sodium	11025610	1
Caelo	Saccharin sodium	12139803	1
Caelo	Saccharin sodium	9381714	1
Caelo	Saccharin sodium	13205414	3
Caelo	Saccharin sodium	14029710	3
Caelo	Saccharin sodium	16221009	4
Caelo	Saccharin sodium	16221010	1
Caelo	Saccharin sodium	17329405	2
Caelo	Saccharin sodium	16221005	11
Caelo	Saccharin sodium	8121801	1
Caelo	Saccharin sodium	14029709	1
Caelo	Saccharin sodium	17329409	2
Caelo	Saccharin sodium	12362607	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Saccharin sodium	18210407	3
Caelo	Saccharin sodium	15087713	16
Caelo	Saccharin sodium	21092015A	1
Caelo	Saccharin sodium	15084101,PZN10220091	1
Caelo	Saccharin sodium	15087704	11
Caelo	Saccharin sodium	15087703	7
Caelo	Saccharin sodium	15084101	5
Caelo	Saccharin sodium	14213405	5
Caelo	Saccharin sodium	13205411	1
Caelo	Saccharin sodium	14213402	2
Caelo	Saccharin sodium	1402970	1
Caelo	Saccharin sodium	14029707	4
Caelo	Saccharin sodium	14029705	6
Caelo	Saccharin sodium	17211002	9
Caelo	Saccharin sodium	17211001	2
Caelo	Saffron filaments	11071310	1
Caelo	Saffron filaments	10184822	1
Caelo	Saffron filaments	12108311	2
Caelo	Saffron filaments	12108308	2
Caelo	Saffron filaments	21011403	1
Caelo	Saffron filaments	10184811	2
Caelo	Saffron filaments	15101305	1
Caelo	Sage leaf, cut	22061201	1
Caelo	Sage leaf, cut	17061302	1
Caelo	Sage leaf, cut	12370011	11
Caelo	Sage leaf, cut	12370009	12
Caelo	Sage leaf, cut	12370006	8
Caelo	Sage leaf, cut	11347513	1
Caelo	Sage leaf, cut	15041304	1
Caelo	Sage leaf, cut	12370003	4
Caelo	Sage leaf, cut	1237003	1
Caelo	Sage leaf, cut	12370001	17
Caelo	Sage leaf, cut	103311513	1
Caelo	Sage leaf, cut	11347509	4
Caelo	Sage leaf, cut	11347511	9
Caelo	Sage leaf, cut	25061302	1
Caelo	Sage leaf, cut	20121201	1
Caelo	Sage leaf, cut	22111203	1
Caelo	Sage leaf, cut	12112101	1
Caelo	Sage leaf, cut	11347508	3
Caelo	Sage leaf, cut	11347503	1
Caelo	Sage leaf, cut	11056125	3
Caelo	Sage leaf, cut	111012004	1
Caelo	Sage leaf, cut	410211513	1
Caelo	Sage leaf, cut	11056123	5
Caelo	Sage leaf, cut	1209211513	1
Caelo	Sage leaf, cut	31081202	1
Caelo	Sage leaf, cut	22081203	1
Caelo	Sage leaf, cut	11056109	5
Caelo	Sage leaf, cut	15121201	1
Caelo	Sage leaf, cut	312010	1
Caelo	Sage leaf, cut	12021302	1
Caelo	Sage leaf, cut	13266606	1
Caelo	Sage leaf, cut	13266626	1
Caelo	Sage leaf, cut	13266621	1
Caelo	Sage leaf, cut	26081301	1
Caelo	Sage leaf, cut	13266605	3
Caelo	Sage leaf, cut	13266603	2
Caelo	Sage leaf, cut	12370014-3541411092	1
Caelo	Sage leaf, cut	12370017	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sage leaf, cut	13011405	1
Caelo	Sage leaf, cut	8011420	1
Caelo	Sage leaf, cut	12370015	2
Caelo	Sage leaf, cut	12370014-3477312161	1
Caelo	Sage leaf, cut	9121302	1
Caelo	Sage leaf, cut	12380009	1
Caelo	Sage leaf, cut	14111315	1
Caelo	Sage leaf, cut	279939	1
Caelo	Sage leaf, cut	12370003-3359319071	1
Caelo	Sage leaf, cut	12370013	3
Caelo	Sage leaf, cut	12370014	8
Caelo	Sage leaf, cut	12370004	2
Caelo	Sage leaf, cut	123770010	1
Caelo	Salicylic acid	160109123	1
Caelo	Salicylic acid	8170717	1
Caelo	Salicylic acid	1601916	1
Caelo	Salicylic acid	16050014	15
Caelo	Salicylic acid	16010922	8
Caelo	Salicylic acid	14076203	8
Caelo	Salicylic acid	13032014A	1
Caelo	Salicylic acid	13344003	11
Caelo	Salicylic acid	7071403	1
Caelo	Salicylic acid	12134719	3
Caelo	Salicylic acid	12134720	8
Caelo	Salicylic acid	12374213	11
Caelo	Salicylic acid	1274210	1
Caelo	Salicylic acid	3601413011-12374211	1
Caelo	Salicylic acid	12374205	9
Caelo	Salicylic acid	1714Q1045	1
Caelo	Salicylic acid	12374212	2
Caelo	Salicylic acid	1334403	1
Caelo	Salicylic acid	14076201	3
Caelo	Salicylic acid	12134721	2
Caelo	Salicylic acid	13344002	20
Caelo	Salicylic acid	16i15-B08-0001	1
Caelo	Salicylic acid	9112013A	1
Caelo	Salicylic acid	16312003	6
Caelo	Salicylic acid	12374210	5
Caelo	Salicylic acid	12374204	2
Caelo	Salicylic acid	160500004	3
Caelo	Salicylic acid	20111306	1
Caelo	Salicylic acid	12134717	9
Caelo	Salicylic acid	12134712	5
Caelo	Salicylic acid	16312005	18
Caelo	Salicylic acid	16072014A	1
Caelo	Salicylic acid	12134710	1
Caelo	Salicylic acid	12374211	18
Caelo	Salicylic acid	123134712	1
Caelo	Salicylic acid	1913Q-01045	1
Caelo	Salicylic acid	10355912	1
Caelo	Salicylic acid	9091301	1
Caelo	Salicylic acid	12134713	1
Caelo	Salicylic acid	12134714	19
Caelo	Salicylic acid	12134715	8
Caelo	Salicylic acid	13B14-No1	1
Caelo	Salicylic acid	1234	1
Caelo	Salicylic acid	1234709	1
Caelo	Salicylic acid	12374207-3454311172	1
Caelo	Salicylic acid	133454002	1
Caelo	Salicylic acid	152117304	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Salicylic acid	14076207	10
Caelo	Salicylic acid	15084401	3
Caelo	Salicylic acid	14300804	8
Caelo	Salicylic acid	16010907	6
Caelo	Salicylic acid	16010903	10
Caelo	Salicylic acid	6/20151706	1
Caelo	Salicylic acid	16050013	17
Caelo	Salicylic acid	1600917	1
Caelo	Salicylic acid	14300809	14
Caelo	Salicylic acid	16010916	2
Caelo	Salicylic acid	14300805	10
Caelo	Salicylic acid	16010906	11
Caelo	Salicylic acid	14336912	8
Caelo	Salicylic acid	23071510	1
Caelo	Salicylic acid	15084403	9
Caelo	Salicylic acid	15084404	4
Caelo	Salicylic acid	151224/1	1
Caelo	Salicylic acid	16010921	1
Caelo	Salicylic acid	15217308	8
Caelo	Salicylic acid	15211808	12
Caelo	Salicylic acid	14336911	20
Caelo	Salicylic acid	3816Q-01045	1
Caelo	Salicylic acid	14076206	6
Caelo	Salicylic acid	16010911	4
Caelo	Salicylic acid	14300806	1
Caelo	Salicylic acid	14076204	11
Caelo	Salicylic acid	13344014	7
Caelo	Salicylic acid	10355907	1
Caelo	Salicylic acid	10112014B	1
Caelo	Salicylic acid	14076211	1
Caelo	Salicylic acid	18121412	1
Caelo	Salicylic acid	14300801	6
Caelo	Salicylic acid	13344013	11
Caelo	Salicylic acid	14076202	7
Caelo	Salicylic acid	13101707	1
Caelo	Salicylic acid	140602701	1
Caelo	Salicylic acid	16010972	1
Caelo	Salicylic acid	16312001	23
Caelo	Salicylic acid	170509AS/16050012	1
Caelo	Salicylic acid	16010918	9
Caelo	Salicylic acid	14300802	6
Caelo	Salicylic acid	1430802	1
Caelo	Salicylic acid	13344018	18
Caelo	Salicylic acid	16010923	6
Caelo	Salicylic acid	13344020	15
Caelo	Salicylic acid	16050016	1
Caelo	Salicylic acid	207102017/16050017	1
Caelo	Salicylic acid	12134702	12
Caelo	Salicylic acid	19061205	1
Caelo	Salicylic acid	15211807	2
Caelo	Salicylic acid	11197123	4
Caelo	Salicylic acid	170301506	1
Caelo	Salicylic acid	12036202	3
Caelo	Salicylic acid	10001612	1
Caelo	Salicylic acid	17123301	1
Caelo	Salicylic acid	18050009	6
Caelo	Salicylic acid	9220715	1
Caelo	Salicylic acid	17051803	1
Caelo	Salicylic acid	11197106	1
Caelo	Salicylic acid	17301506	23

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Salicylic acid	16010901	3
Caelo	Salicylic acid	16312008	1
Caelo	Salicylic acid	17133307	23
Caelo	Salicylic acid	102455-0001	1
Caelo	Salicylic acid	1713301	1
Caelo	Salicylic acid	17301504	2
Caelo	Salicylic acid	11311922	13
Caelo	Salicylic acid	17301503	11
Caelo	Salicylic acid	17301512	7
Caelo	Salicylic acid	17301502	9
Caelo	Salicylic acid	17301511	1
Caelo	Salicylic acid	18050002	7
Caelo	Salicylic acid	1310516	1
Caelo	Salicylic acid	16010904	9
Caelo	Salicylic acid	10355916	1
Caelo	Salicylic acid	10291803	1
Caelo	Salicylic acid	15081601	1
Caelo	Salicylic acid	11311902	2
Caelo	Salicylic acid	11311903	1
Caelo	Salicylic acid	14336920	21
Caelo	Salicylic acid	16050003	2
Caelo	Salicylic acid	11197117	1
Caelo	Salicylic acid	11311906	3
Caelo	Salicylic acid	14336921	1
Caelo	Salicylic acid	14336919	19
Caelo	Salicylic acid	11311910	2
Caelo	Salicylic acid	17301509	14
Caelo	Salicylic acid	18050007	8
Caelo	Salicylic acid	16050006	1
Caelo	Salicylic acid	18F05-B02-352210	3
Caelo	Salicylic acid	180104201	5
Caelo	Salicylic acid	18050006	11
Caelo	Salicylic acid	4515E-01045	1
Caelo	Salicylic acid	16010905	10
Caelo	Salicylic acid	17133303	15
Caelo	Salicylic acid	18147302	2
Caelo	Salicylic acid	20270312-2	1
Caelo	Salicylic acid	163120	1
Caelo	Salicylic acid	12036217	2
Caelo	Salicylic acid	12036223	8
Caelo	Salicylic acid	12134709	22
Caelo	Salicylic acid	10355911	2
Caelo	Salicylic acid	25041304	1
Caelo	Salicylic acid	4010Q-01045	1
Caelo	Salicylic acid	16050004	12
Caelo	Salicylic acid	11311923	2
Caelo	Salicylic acid	12036315	1
Caelo	Salicylic acid	12323402	2
Caelo	Salicylic acid	12036225	4
Caelo	Salicylic acid	18147307	1
Caelo	Salicylic acid	17133301	8
Caelo	Salicylic acid	190263	1
Caelo	Salicylic acid	12323401	2
Caelo	Salicylic acid	72(18230009)	1
Caelo	Salicylic acid	18230002	1
Caelo	Salicylic acid	12134711	9
Caelo	Salicylic acid	16312005, 20171127Ma4	1
Caelo	Salicylic acid	A-100513-1	1
Caelo	Salicylic acid	11311905	1
Caelo	Salicylic acid	16312010	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Salicylic acid	15217301	9
Caelo	Salicylic acid	160405AS	1
Caelo	Salicylic acid	18147303	3
Caelo	Salicylic acid	20151112-4	1
Caelo	Salicylic acid	17301612	1
Caelo	Salicylic acid	11197118	1
Caelo	Salicylic acid	18230004	4
Caelo	Salicylic acid	18147310	2
Caelo	Salicylic acid	1197106	1
Caelo	Salicylic acid		3
Caelo	Salicylic acid	12036211	4
Caelo	Salicylic acid	7466	1
Caelo	Salicylic acid	12036209	4
Caelo	Salicylic acid	12036210	3
Caelo	Salicylic acid	17133304	1
Caelo	Salicylic acid	17301505	14
Caelo	Salicylic acid	16312004	1
Caelo	Salicylic acid	12036208	1
Caelo	Salicylic acid	3111J01045	1
Caelo	Salicylic acid	18230007	6
Caelo	Salicylic acid	16312009	7
Caelo	Salicylic acid	1738980	1
Caelo	Salicylic acid	14336907	23
Caelo	Salicylic acid	11197103	1
Caelo	Salicylic acid	14336908	2
Caelo	Salicylic acid	15217305	6
Caelo	Salicylic acid	15211810	9
Caelo	Salicylic acid	Val.Mes.:2315Q01045	1
Caelo	Salicylic acid	1045	1
Caelo	Salicylic acid	16010917	15
Caelo	Salicylic acid	14336904	2
Caelo	Salicylic acid	15211803	1
Caelo	Salicylic acid	16010913	13
Caelo	Salicylic acid	16050005	11
Caelo	Salicylic acid	16010908	3
Caelo	Salicylic acid	16E02-B09-325355	1
Caelo	Salicylic acid	16010912	4
Caelo	Salicylic acid	16050011	11
Caelo	Salicylic acid	15217302	14
Caelo	Salicylic acid	15211809	12
Caelo	Salicylic acid	Validierungsmessung	1
Caelo	Salicylic acid	15211806	3
Caelo	Salicylic acid	15217303	9
Caelo	Salicylic acid	14336906	3
Caelo	Salicylic acid	15217304	4
Caelo	Salicylic acid	14336905	13
Caelo	Sarsaparilla root, cut	13079201	1
Caelo	Savory, rubbed	12200801	1
Caelo	Savory, rubbed	13172404	2
Caelo	Savory, rubbed	12177203	7
Caelo	Savory, rubbed	12177201	1
Caelo	Savory, rubbed	12177207	2
Caelo	Selenium yeast	15276208	1
Caelo	Selenium yeast	13347906	1
Caelo	Selenium yeast	13347911	1
Caelo	Senna leaf, cut	12115501	5
Caelo	Senna leaf, cut	2909211514	1
Caelo	Senna leaf, cut	12273106	3
Caelo	Senna leaf, cut	12273104	2
Caelo	Senna leaf, cut	12273107	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Senna leaf, cut	1273103	1
Caelo	Senna leaf, cut	2331B120521	4
Caelo	Senna leaf, cut	12273111	10
Caelo	Senna leaf, cut	711211514	1
Caelo	Senna leaf, cut	10126502	1
Caelo	Senna leaf, cut	131316701	1
Caelo	Senna leaf, cut	12115503	6
Caelo	Senna leaf, cut	15101303	1
Caelo	Senna leaf, cut	103311514	2
Caelo	Senna leaf, cut	11041009	1
Caelo	Senna leaf, cut	12273105	6
Caelo	Senna leaf, cut	12273114	2
Caelo	Senna leaf, cut	12273113	3
Caelo	Senna leaf, cut	2231B13014	1
Caelo	Senna leaf, cut	2413Q01514	1
Caelo	Senna leaf, cut	12273119	4
Caelo	Senna leaf, cut	12115502	1
Caelo	Senna leaf, cut	12273101	4
Caelo	Senna leaf, cut	9071302	1
Caelo	Senna leaf, cut	12273115	11
Caelo	Senna leaf, cut	1803311514	1
Caelo	Senna leaf, cut	12273120	3
Caelo	Senna leaf, cut	12273118	3
Caelo	Senna leaf, cut	313082	1
Caelo	Senna leaf, cut	5091306	1
Caelo	Senna leaf, cut	3813E1514	1
Caelo	Senna leaf, cut	13136703	3
Caelo	Senna leaf, cut	13136704	13
Caelo	Senna leaf, cut	1316703	2
Caelo	Senna leaf, cut	4113Q1514	1
Caelo	Senna leaf, cut	11277301	1
Caelo	Senna leaf, cut	13136705	3
Caelo	Senna leaf, cut	111113a	1
Caelo	Senna leaf, cut	12276111	1
Caelo	Senna leaf, cut	13136706	1
Caelo	Senna leaf, cut	1316706	1
Caelo	Senna leaf, cut	13136705-3489312103	1
Caelo	Senna leaf, cut	18071201	1
Caelo	Senna leaf, cut	14306310	1
Caelo	Senna leaf, cut	806211514	1
Caelo	Senna leaf, cut	11277304	4
Caelo	Senna leaf, cut	5013U01514	1
Caelo	Senna leaf, cut	14183805	2
Caelo	Senna leaf, cut	13136724	1
Caelo	Senna leaf, cut	12369903	1
Caelo	Senna leaf, cut	514R1514	1
Caelo	Senna leaf, cut	21011401	1
Caelo	Senna leaf, cut	11277305	1
Caelo	Senna leaf, cut	14306324	1
Caelo	Senna leaf, cut	11277307	1
Caelo	Senna leaf, cut	13136718	1
Caelo	Senna leaf, cut	11277316	4
Caelo	Senna leaf, cut	14306311	2
Caelo	Senna leaf, powdered	11193307	1
Caelo	Senna leaf, powdered	17223202	1
Caelo	Senna leaf, powdered	14347606	1
Caelo	Senna pods, tinnevelly	11003013	1
Caelo	Senna pods, tinnevelly	12267602	2
Caelo	Senna pods, tinnevelly	12267201	1
Caelo	Senna pods, tinnevelly	11003004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Senna pods, tinnevelly	110030009	1
Caelo	Senna pods, tinnevelly	12267610	1
Caelo	Senna pods, tinnevelly	12267608	1
Caelo	Senna pods, tinnevelly	12267609	3
Caelo	Senna pods, tinnevelly	12267007	1
Caelo	Senna pods, tinnevelly	12267604	1
Caelo	Senna pods, tinnevelly	12267605	1
Caelo	Senna pods, tinnevelly	122676707	1
Caelo	Senna pods, tinnevelly	11003009	2
Caelo	Senna pods, tinnevelly	12267606	1
Caelo	Senna pods, tinnevelly	22061202	1
Caelo	Shepherd's purse, cut	11036410	2
Caelo	Shepherd's purse, cut	14061201	1
Caelo	Shepherd's purse, cut	11036401	1
Caelo	Shepherd's purse, cut	11091202	1
Caelo	Shepherd's purse, cut	11328205	1
Caelo	Shepherd's purse, cut	11328209	6
Caelo	Shepherd's purse, cut	12354201	11
Caelo	Shepherd's purse, cut	15101202	1
Caelo	Shepherd's purse, cut	269595	1
Caelo	Shepherd's purse, cut	11328204	3
Caelo	Shepherd's purse, cut	11328208	6
Caelo	Shepherd's purse, cut	11328207	8
Caelo	Shepherd's purse, cut	11328203	5
Caelo	Shepherd's purse, cut	11328201	6
Caelo	Shepherd's purse, cut	11328206	2
Caelo	Shepherd's purse, cut	10061301	1
Caelo	Shepherd's purse, cut	11036404	1
Caelo	Shepherd's purse, cut	11036406	5
Caelo	Shepherd's purse, cut	12354203-3396310111	1
Caelo	Shepherd's purse, cut	12354205	6
Caelo	Shepherd's purse, cut	12354001	1
Caelo	Shepherd's purse, cut	3121301	1
Caelo	Shepherd's purse, cut	12254202	1
Caelo	Shepherd's purse, cut	12354203	13
Caelo	Shepherd's purse, cut	12354206	4
Caelo	Shepherd's purse, cut	12354207	2
Caelo	Shepherd's purse, cut	12354204	16
Caelo	Shepherd's purse, cut	12354209	1
Caelo	Shepherd's purse, cut	12354202	2
Caelo	Shepherd's purse, cut	1132820111-14	1
Caelo	Shepherd's purse, cut	13344203	1
Caelo	Shepherd's purse, cut	12354503	1
Caelo	Sildenafil citrate	15097303	1
Caelo	Sildenafil citrate	101051528	1
Caelo	Sildenafil citrate	16049819	3
Caelo	Sildenafil citrate	16049808	5
Caelo	Sildenafil citrate	16049814	9
Caelo	Sildenafil citrate	2101403	1
Caelo	Sildenafil citrate	16049812	2
Caelo	Sildenafil citrate	14074709	1
Caelo	Sildenafil citrate	15097312	2
Caelo	Sildenafil citrate	14074707	2
Caelo	Sildenafil citrate	15097311	1
Caelo	Sildenafil citrate	15097310	1
Caelo	Sildenafil citrate	14074703	2
Caelo	Sildenafil citrate	15097317	1
Caelo	Sildenafil citrate	16049802	10
Caelo	Sildenafil citrate	16049821	8
Caelo	Sildenafil citrate	16049820	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sildenafil citrate	15097307	2
Caelo	Sildenafil citrate	17071504	7
Caelo	Sildenafil citrate	16049826	1
Caelo	Sildenafil citrate	17071502	2
Caelo	Sildenafil citrate	15F11-B02-309081	1
Caelo	Sildenafil citrate	15F11-B02-314725	2
Caelo	Sildenafil citrate	15097301	2
Caelo	Sildenafil citrate	16049823	2
Caelo	Sildenafil citrate	17344008	1
Caelo	Sildenafil citrate	14074715	5
Caelo	Sildenafil citrate	17344006	1
Caelo	Sildenafil citrate	16049824	1
Caelo	Sildenafil citrate	1407474	1
Caelo	Sildenafil citrate	17344001	1
Caelo	Sildenafil citrate	17344014	1
Caelo	Sildenafil citrate	14074714	5
Caelo	Sildenafil citrate	16049801	1
Caelo	Sildenafil citrate	14074716	1
Caelo	Sildenafil citrate	15097315	1
Caelo	Sildenafil citrate	16049810	2
Caelo	Sildenafil citrate	15097306	1
Caelo	Silverweed, cut	11328110	1
Caelo	Silverweed, cut	11328112	2
Caelo	Silverweed, cut	11016404	1
Caelo	Silverweed, cut	20081202	1
Caelo	Silverweed, cut	11016407	1
Caelo	Silverweed, cut	8111204	1
Caelo	Silverweed, cut	11328105	7
Caelo	Silverweed, cut	11328109	13
Caelo	Sodium acetate trihyd...	15145204	6
Caelo	Sodium acetate trihyd...	12027202	1
Caelo	Sodium acetate trihyd...	18006801	5
Caelo	Sodium acetate trihyd...	13411418	1
Caelo	Sodium acetate trihyd...	16300711	7
Caelo	Sodium acetate trihyd...	12027217	1
Caelo	Sodium acetate trihyd...	13295606	1
Caelo	Sodium acetate trihyd...	16300706	5
Caelo	Sodium acetate trihyd...	10353017	1
Caelo	Sodium acetate trihyd...	12316501	1
Caelo	Sodium acetate trihyd...	15118701	3
Caelo	Sodium acetate trihyd...	13411418	1
Caelo	Sodium acetate trihyd...	15145209	6
Caelo	Sodium acetate trihyd...	16300712	1
Caelo	Sodium acetate trihyd...	151452009	1
Caelo	Sodium acetate trihyd...	15145218	2
Caelo	Sodium acetate trihyd...	12027228	1
Caelo	Sodium acetate trihyd...	15145205	1
Caelo	Sodium ascorbate	27091302	1
Caelo	Sodium ascorbate	14079404	1
Caelo	Sodium ascorbate	15180505	1
Caelo	Sodium ascorbate	1604001-01	1
Caelo	Sodium ascorbate	16285710	1
Caelo	Sodium benzoate	17055801	5
Caelo	Sodium benzoate	15044102	5
Caelo	Sodium benzoate	14328702	1
Caelo	Sodium benzoate	16137112	6
Caelo	Sodium benzoate	91443209	2
Caelo	Sodium benzoate	161371	1
Caelo	Sodium benzoate	1503022-01	1
Caelo	Sodium benzoate	16137104	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium benzoate	16137118	7
Caelo	Sodium benzoate	11267016	1
Caelo	Sodium benzoate	1610037-01	3
Caelo	Sodium benzoate	4315A-01891	1
Caelo	Sodium benzoate	16137103	9
Caelo	Sodium benzoate	18055805	1
Caelo	Sodium benzoate	13178647	2
Caelo	Sodium benzoate	17055805	7
Caelo	Sodium benzoate	19091801	1
Caelo	Sodium benzoate	18161201	2
Caelo	Sodium benzoate	17353006	2
Caelo	Sodium benzoate	19054104	1
Caelo	Sodium benzoate	1209018-0	1
Caelo	Sodium benzoate	12158401	1
Caelo	Sodium benzoate	12158419	4
Caelo	Sodium benzoate	15363401	10
Caelo	Sodium benzoate	12158424	1
Caelo	Sodium benzoate	17061306	1
Caelo	Sodium benzoate	15135001	6
Caelo	Sodium benzoate	12158418	1
Caelo	Sodium benzoate	12158439	3
Caelo	Sodium benzoate	13178603	1
Caelo	Sodium benzoate	13178605	3
Caelo	Sodium benzoate	5211I-01891	1
Caelo	Sodium benzoate	13178613	1
Caelo	Sodium benzoate	15044104	1
Caelo	Sodium benzoate	13178619	2
Caelo	Sodium benzoate	13178603	1
Caelo	Sodium benzoate	13178635	6
Caelo	Sodium benzoate	13178627	1
Caelo	Sodium benzoate	1506019-01	6
Caelo	Sodium benzoate	17353002	1
Caelo	Sodium benzoate	15363409	2
Caelo	Sodium benzoate	17055806	1
Caelo	Sodium benzoate	24021501	1
Caelo	Sodium benzoate	10544104	1
Caelo	Sodium benzoate	5117A-01891	1
Caelo	Sodium benzoate	17137604	3
Caelo	Sodium benzoate	15363413	3
Caelo	Sodium benzoate	17137606	15
Caelo	Sodium benzoate	151208NB	1
Caelo	Sodium benzoate	14166714	1
Caelo	Sodium benzoate	18161211	3
Caelo	Sodium benzoate	15363424	10
Caelo	Sodium benzoate	17353001	4
Caelo	Sodium benzoate	2215Q-01891	1
Caelo	Sodium benzoate	15363414	1
Caelo	Sodium benzoate	NB20180313	1
Caelo	Sodium benzoate	17137609	1
Caelo	Sodium benzoate	11346109	3
Caelo	Sodium benzoate	15044104	7
Caelo	Sodium carbonate deca...	11254706	1
Caelo	Sodium carbonate deca...	12239806	1
Caelo	Sodium carbonate deca...	12239810	2
Caelo	Sodium carbonate deca...	80136179	1
Caelo	Sodium carbonate deca...	13036902	2
Caelo	Sodium carbonate deca...	13036907	1
Caelo	Sodium carbonate deca...	13036909	1
Caelo	Sodium carbonate deca...	13036910	1
Caelo	Sodium carbonate deca...	11254708	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium carbonate deca...	80136428	1
Caelo	Sodium citrate	13089329	1
Caelo	Sodium citrate	1234	3
Caelo	Sodium citrate	30081802	1
Caelo	Sodium citrate	4111Q-01900	1
Caelo	Sodium citrate	17323608	15
Caelo	Sodium citrate	171006-02	1
Caelo	Sodium citrate	1123588946	1
Caelo	Sodium citrate	123456789	1
Caelo	Sodium citrate	2413A-01900	1
Caelo	Sodium citrate	140510/1	1
Caelo	Sodium citrate	17323609	1
Caelo	Sodium citrate	12045002	1
Caelo	Sodium citrate	14369424	4
Caelo	Sodium citrate	123tst	1
Caelo	Sodium citrate	11304311	3
Caelo	Sodium citrate	2543	1
Caelo	Sodium citrate	11219012	1
Caelo	Sodium citrate	14369415	1
Caelo	Sodium citrate	14112518	1
Caelo	Sodium citrate	Test	2
Caelo	Sodium citrate	17081501	1
Caelo	Sodium citrate	25383143	1
Caelo	Sodium citrate	11081514	1
Caelo	Sodium citrate	14369407	2
Caelo	Sodium citrate	50815PMS01	1
Caelo	Sodium citrate	14369407	32
Caelo	Sodium citrate	14369413	14
Caelo	Sodium citrate	1234Test	3
Caelo	Sodium citrate	14369411	3
Caelo	Sodium citrate	14112502	7
Caelo	Sodium citrate	1223test	1
Caelo	Sodium citrate	14369416	1
Caelo	Sodium citrate	123456	1
Caelo	Sodium citrate	14369423	11
Caelo	Sodium citrate	123658	1
Caelo	Sodium citrate	2345678	1
Caelo	Sodium citrate	140731NC	1
Caelo	Sodium citrate	16095702	3
Caelo	Sodium citrate	2815A-01900	1
Caelo	Sodium citrate	15212705	4
Caelo	Sodium citrate	14369423	1
Caelo	Sodium citrate	19121802	1
Caelo	Sodium citrate	12045019	5
Caelo	Sodium citrate	TestmessungnachRep	1
Caelo	Sodium citrate	1151351	2
Caelo	Sodium citrate	12045026	10
Caelo	Sodium citrate	1732608	1
Caelo	Sodium citrate	125874566	1
Caelo	Sodium citrate	154662	1
Caelo	Sodium citrate	12045024	2
Caelo	Sodium citrate	128896	1
Caelo	Sodium citrate	18208302	2
Caelo	Sodium citrate	18208307	4
Caelo	Sodium citrate	14369414	1
Caelo	Sodium citrate	16217901	1
Caelo	Sodium citrate	12045024	4
Caelo	Sodium citrate	1234	1
Caelo	Sodium citrate	14838	1
Caelo	Sodium citrate	K-280313-1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium citrate	16217902	1
Caelo	Sodium citrate	1256825	1
Caelo	Sodium citrate	17323612	1
Caelo	Sodium citrate	459842	1
Caelo	Sodium citrate	16095707	6
Caelo	Sodium citrate	21121203	1
Caelo	Sodium citrate	160704001	1
Caelo	Sodium citrate	123test	4
Caelo	Sodium citrate	1732361	1
Caelo	Sodium citrate	11121303	1
Caelo	Sodium citrate	12045016-3451311172	1
Caelo	Sodium citrate	49131-01900	1
Caelo	Sodium citrate	17323601	18
Caelo	Sodium citrate	13/24-B13-287921	2
Caelo	Sodium citrate	14369406	1
Caelo	Sodium citrate	11219012	1
Caelo	Sodium citrate	8/AKD41	1
Caelo	Sodium citrate	16217913	1
Caelo	Sodium citrate	70214PJK01	1
Caelo	Sodium citrate	16217912	3
Caelo	Sodium citrate	16217914	3
Caelo	Sodium citrate	13089319	3
Caelo	Sodium citrate	1411502	1
Caelo	Sodium citrate	14112502	12
Caelo	Sodium citrate	16217908	1
Caelo	Sodium citrate	131103501	1
Caelo	Sodium citrate	13089335	7
Caelo	Sodium citrate	123Test	1
Caelo	Sodium citrate	12045001	1
Caelo	Sodium citrate	12045008	3
Caelo	Sodium citrate	13089331	7
Caelo	Sodium citrate	16217906	1
Caelo	Sodium citrate	12897	1
Caelo	Sodium citrate	171025NC/16217909	1
Caelo	Sodium citrate	13089305	1
Caelo	Sodium citrate	13089317	7
Caelo	Sodium citrate	13089325	3
Caelo	Sodium citrate	1308908	1
Caelo	Sodium citrate	17H7-B04-345098	1
Caelo	Sodium citrate	12045016	5
Caelo	Sodium citrate	4680	1
Caelo	Sodium citrate	16217905	13
Caelo	Sodium citrate	16F10-B01-326093	1
Caelo	Sodium citrate	12045012	1
Caelo	Sodium citrate	220914PML01	1
Caelo	Sodium citrate	11219012	1
Caelo	Sodium citrate	16217909	6
Caelo	Sodium citrate	14360407	1
Caelo	Sodium citrate	13089314	1
Caelo	Sodium citrate	14369402	2
Caelo	Sodium citrate	12045006	1
Caelo	Sodium citrate	12045011	6
Caelo	Sodium citrate	13089319	1
Caelo	Sodium citrate	14112520	1
Caelo	Sodium citrate	14112517	4
Caelo	Sodium citrate	123Test	5
Caelo	Sodium citrate	1042	1
Caelo	Sodium citrate	13089308	8
Caelo	Sodium citrate	14112502	5
Caelo	Sodium citrate	12045004	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium citrate	1161207	1
Caelo	Sodium citrate	14112517	3
Caelo	Sodium citrate	14112517	2
Caelo	Sodium citrate	13089304	3
Caelo	Sodium cromoglicate	15224404	1
Caelo	Sodium cromoglicate	5113A-02831	1
Caelo	Sodium cromoglicate	16270705	2
Caelo	Sodium cromoglicate	12301102	1
Caelo	Sodium cyclamate	15201606	2
Caelo	Sodium cyclamate	11204205	1
Caelo	Sodium cyclamate	16119310	1
Caelo	Sodium cyclamate	17041102	5
Caelo	Sodium cyclamate	17362202	4
Caelo	Sodium cyclamate	15201614	1
Caelo	Sodium cyclamate	16119308	3
Caelo	Sodium cyclamate	18216702	2
Caelo	Sodium cyclamate	15201608	2
Caelo	Sodium cyclamate	16119302	5
Caelo	Sodium cyclamate	16119302/2280217	1
Caelo	Sodium cyclamate	179617362207	1
Caelo	Sodium cyclamate	10380915	2
Caelo	Sodium cyclamate	15201601	2
Caelo	Sodium dihydrogen pho...	12198517	2
Caelo	Sodium dihydrogen pho...	13D08-N25	1
Caelo	Sodium dihydrogen pho...	12198506	1
Caelo	Sodium dihydrogen pho...	12198510	2
Caelo	Sodium dodecyl sulfate	12153010	3
Caelo	Sodium dodecyl sulfate	12153019	1
Caelo	Sodium dodecyl sulfate	13014604	2
Caelo	Sodium dodecyl sulfate	13099006	3
Caelo	Sodium dodecyl sulfate	13099011	3
Caelo	Sodium dodecyl sulfate	13411302	2
Caelo	Sodium dodecyl sulfate	13411305	3
Caelo	Sodium dodecyl sulfate	13411304	1
Caelo	Sodium dodecyl sulfate	14127406	1
Caelo	Sodium dodecyl sulfate	14127411	1
Caelo	Sodium dodecyl sulfate	14283904	3
Caelo	Sodium dodecyl sulfate	18045202	1
Caelo	Sodium dodecyl sulfate	11220208	1
Caelo	Sodium dodecyl sulfate	17095511	4
Caelo	Sodium dodecyl sulfate	16229220	1
Caelo	Sodium dodecyl sulfate	17095504	3
Caelo	Sodium dodecyl sulfate	16229218	1
Caelo	Sodium dodecyl sulfate	16229209	1
Caelo	Sodium dodecyl sulfate	16229217	4
Caelo	Sodium dodecyl sulfate	16229213	2
Caelo	Sodium dodecyl sulfate	16229206	4
Caelo	Sodium dodecyl sulfate	15344515	2
Caelo	Sodium dodecyl sulfate	15344511	1
Caelo	Sodium dodecyl sulfate	15344505	2
Caelo	Sodium dodecyl sulfate	15344503	3
Caelo	Sodium dodecyl sulfate	15119907	3
Caelo	Sodium dodecyl sulfate	15095404	1
Caelo	Sodium dodecyl sulfate	12153004	1
Caelo	Sodium dodecyl sulfate	12153003	1
Caelo	Sodium dodecyl sulfate	26031508	1
Caelo	Sodium dodecyl sulfate	14224206	1
Caelo	Sodium edetate	15116602	4
Caelo	Sodium edetate	15035603	2
Caelo	Sodium edetate	14270506	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium edetate	92159379	1
Caelo	Sodium edetate	13142002	1
Caelo	Sodium edetate	12181704	1
Caelo	Sodium edetate	1427006	1
Caelo	Sodium edetate	14270505	1
Caelo	Sodium edetate	15427714	20
Caelo	Sodium edetate	13E15-N05	1
Caelo	Sodium edetate	12181706	2
Caelo	Sodium edetate	14270505	2
Caelo	Sodium edetate	15427717	8
Caelo	Sodium edetate	1517E-01925	1
Caelo	Sodium edetate	18007110	4
Caelo	Sodium edetate	10887102	1
Caelo	Sodium edetate	18007102	3
Caelo	Sodium edetate	16136404	5
Caelo	Sodium edetate	16337302	6
Caelo	Sodium edetate	16337316	5
Caelo	Sodium edetate	3211I-01925	1
Caelo	Sodium edetate	17336102	9
Caelo	Sodium edetate	12181703	3
Caelo	Sodium edetate	14031801	1
Caelo	Sodium edetate	11074703	1
Caelo	Sodium edetate	12181711	2
Caelo	Sodium edetate	15116603	14
Caelo	Sodium edetate	13395007	2
Caelo	Sodium edetate	13/23-N03-3472312131	1
Caelo	Sodium edetate	14270505,PZN01397301	1
Caelo	Sodium edetate	15427703	15
Caelo	Sodium edetate	13142014	1
Caelo	Sodium edetate	11074608	2
Caelo	Sodium edetate	13142010	2
Caelo	Sodium edetate	13142015	3
Caelo	Sodium edetate	12181702	1
Caelo	Sodium edetate	2101414	1
Caelo	Sodium edetate	13142009	4
Caelo	Sodium edetate	15035604	4
Caelo	Sodium edetate	13142002	3
Caelo	Sodium edetate	15427705	5
Caelo	Sodium edetate	13142004	5
Caelo	Sodium edetate	2913M1925	1
Caelo	Sodium edetate	13142015	1
Caelo	Sodium edetate	12254507	1
Caelo	Sodium edetate	13395002	3
Caelo	Sodium edetate	3913A-01925	1
Caelo	Sodium edetate	13E15-N05	1
Caelo	Sodium edetate	15427710	6
Caelo	Sodium edetate	12181714	3
Caelo	Sodium edetate	12181714	2
Caelo	Sodium edetate	13142008	3
Caelo	Sodium edetate	13395007	2
Caelo	Sodium edetate	13142008	4
Caelo	Sodium edetate	17d25-b08-339881	1
Caelo	Sodium hydrogen carbo...	13234202	3
Caelo	Sodium hydrogen carbo...	13234007	4
Caelo	Sodium hydrogen carbo...	12324708	1
Caelo	Sodium hydrogen carbo...	8101203	1
Caelo	Sodium hydrogen carbo...	12132203	3
Caelo	Sodium hydrogen carbo...	12324725	1
Caelo	Sodium hydrogen carbo...	13060103	8
Caelo	Sodium hydrogen carbo...	12356303	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium hydrogen carbo...	1203030-1	1
Caelo	Sodium hydrogen carbo...	3012Q-01892	1
Caelo	Sodium hydrogen carbo...	12109204	3
Caelo	Sodium hydrogen carbo...	11304802	2
Caelo	Sodium hydrogen carbo...	11304809	1
Caelo	Sodium hydrogen carbo...	12132219	8
Caelo	Sodium hydrogen carbo...	12132204	7
Caelo	Sodium hydrogen carbo...	1219204	1
Caelo	Sodium hydrogen carbo...	12324702	7
Caelo	Sodium hydrogen carbo...	1234702	2
Caelo	Sodium hydrogen carbo...	11304813	1
Caelo	Sodium hydrogen carbo...	13060117	1
Caelo	Sodium hydrogen carbo...	12356401	3
Caelo	Sodium hydrogen carbo...	13060120	5
Caelo	Sodium hydrogen carbo...	13390520	1
Caelo	Sodium hydrogen carbo...	13186721	1
Caelo	Sodium hydrogen carbo...	13234005	4
Caelo	Sodium hydrogen carbo...	13234205	1
Caelo	Sodium hydrogen carbo...	12132503	1
Caelo	Sodium hydrogen carbo...	13390505	1
Caelo	Sodium hydrogen carbo...	6011404	1
Caelo	Sodium hydrogen carbo...	13186711	1
Caelo	Sodium hydrogen carbo...	13234004	2
Caelo	Sodium hydrogen carbo...	15328602	1
Caelo	Sodium hydrogen carbo...	12356407	2
Caelo	Sodium hydrogen carbo...	12356403	1
Caelo	Sodium hydrogen carbo...	13186704	1
Caelo	Sodium hydrogen carbo...	13101706	1
Caelo	Sodium hydrogen carbo...	12109205	6
Caelo	Sodium hydrogen carbo...	13186706	1
Caelo	Sodium hydroxide pell...	11070821	1
Caelo	Sodium hydroxide pell...	12368105	1
Caelo	Sodium hydroxide pell...	15312102	1
Caelo	Sodium methyl-4-hydro...	17051903	1
Caelo	Sodium methyl-4-hydro...	18060801	1
Caelo	Sodium monohydrogenph...	12345025	1
Caelo	Sodium monohydrogenph...	71199028	1
Caelo	Sodium monohydrogenph...	1169315	1
Caelo	Sodium monohydrogenph...	12345016	7
Caelo	Sodium monohydrogenph...	AKD7	1
Caelo	Sodium monohydrogenph...	9302613	2
Caelo	Sodium monohydrogenph...	11169315	1
Caelo	Sodium monohydrogenph...	11169303	2
Caelo	Sodium monohydrogenph...	11169309	3
Caelo	Sodium monohydrogenph...	12345001	2
Caelo	Sodium monohydrogenph...	12256602	2
Caelo	Sodium monohydrogenph...	12345013	6
Caelo	Sodium monohydrogenph...	12345004	7
Caelo	Sodium monohydrogenph...	83002059	1
Caelo	Sodium monohydrogenph...	12345007	3
Caelo	Sodium monohydrogenph...	13028103	1
Caelo	Sodium monohydrogenph...	72006059	1
Caelo	Sodium monohydrogenph...	13028003	1
Caelo	Sodium monohydrogenph...	12109109	1
Caelo	Sodium monohydrogenph...	12109002	2
Caelo	Sodium monohydrogenph...	12109106	1
Caelo	Sodium monohydrogenph...	12109108	2
Caelo	Sodium monohydrogenph...	11109513	1
Caelo	Sodium monohydrogenph...	13028001	1
Caelo	Sodium monohydrogenph...	10135304	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sodium perborate, hyd...	13050807	1
Caelo	Sodium perborate, hyd...	13050805	1
Caelo	Sodium perborate, hyd...	13050801	2
Caelo	Sodium salicylate, po...	18061305	1
Caelo	Sodium sulfate decahy...	12157001	1
Caelo	Sodium sulfate decahy...	18101202	1
Caelo	Sodium sulfate decahy...	11012207	1
Caelo	Sodium sulfate decahy...	11012210	2
Caelo	Sodium sulfate decahy...	12256902	1
Caelo	Sodium sulfate decahy...	12257901	1
Caelo	Sodium sulfate decahy...	12157013	1
Caelo	Sodium sulfate decahy...	12157007	1
Caelo	Sodium sulfate decahy...	12157006	3
Caelo	Sodium sulfate decahy...	70706387	1
Caelo	Sodium sulfate decahy...	205026	1
Caelo	Sodium sulfate decahy...	12157012	18
Caelo	Sodium sulfate decahy...	12157011	1
Caelo	Sodium sulfate decahy...	61599506	1
Caelo	Sodium sulfate decahy...	12157005	13
Caelo	Sodium tetraborate	16081801	1
Caelo	Sodium tetraborate	13290118	1
Caelo	Sodium tetraborate	11369431	1
Caelo	Sodium tetraborate	14310905	2
Caelo	Sodium tetraborate	16137201	1
Caelo	Sodium tetraborate	21051402	1
Caelo	Sodium thiosulfate pe...	26041302	1
Caelo	Sodium thiosulfate pe...	80426288	1
Caelo	Softisan® 601	13296104	1
Caelo	Softisan® 601	14099608	1
Caelo	Softisan® 601	15184901	1
Caelo	Solanum dulcamara ste...	13254301	1
Caelo	Solanum dulcamara ste...	132564301	1
Caelo	Solanum dulcamara ste...	10375605	1
Caelo	Solanum dulcamara ste...	103745604	1
Caelo	Sorbic acid	17343422	7
Caelo	Sorbic acid	21041801	1
Caelo	Sorbic acid	18287605	4
Caelo	Sorbic acid	10305209	1
Caelo	Sorbic acid	1827669	1
Caelo	Sorbic acid	162508016	1
Caelo	Sorbic acid	18284605	1
Caelo	Sorbic acid	15123912	1
Caelo	Sorbic acid	18287609	1
Caelo	Sorbic acid	17343411	13
Caelo	Sorbic acid	5012E-01048	1
Caelo	Sorbic acid	11203806	2
Caelo	Sorbic acid	11203807	10
Caelo	Sorbic acid	29071201	1
Caelo	Sorbic acid	16208016	5
Caelo	Sorbic acid	16208009	11
Caelo	Sorbic acid	28041726	1
Caelo	Sorbic acid	170427AS/16208005	1
Caelo	Sorbic acid	1620808005	1
Caelo	Sorbic acid	16208005	24
Caelo	Sorbic acid	15382010	14
Caelo	Sorbic acid	15382001	1
Caelo	Sorbic acid	14380715	1
Caelo	Sorbic acid	15382002	10
Caelo	Sorbic acid	14176209	1
Caelo	Sorbic acid	14380714	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sorbic acid	14380706	1
Caelo	Sorbic acid	14380713	2
Caelo	Sorbic acid	12330608	1
Caelo	Sorbic acid	14380714	7
Caelo	Sorbic acid	14380706	6
Caelo	Sorbic acid	14176212	3
Caelo	Sorbic acid	14380704	1
Caelo	Sorbic acid	14380702	1
Caelo	Sorbic acid	14176214	6
Caelo	Sorbic acid	14176213	1
Caelo	Sorbic acid	14176214	3
Caelo	Sorbic acid	14380702	2
Caelo	Sorbic acid	14176212	5
Caelo	Sorbic acid	14176207	1
Caelo	Sorbic acid	14176209	3
Caelo	Sorbic acid	13334010	7
Caelo	Sorbic acid	12038009	1
Caelo	Sorbic acid	13334003	4
Caelo	Sorbic acid	4811Q-01048	1
Caelo	Sorbic acid	16208012	6
Caelo	Sorbic acid	16208011	9
Caelo	Sorbic acid	17343401	4
Caelo	Sorbic acid	16208020	7
Caelo	Sorbic acid	10305208	1
Caelo	Sorbic acid	11203812	1
Caelo	Sorbic acid	12330618	13
Caelo	Sorbic acid	162080	1
Caelo	Sorbic acid	12330608	8
Caelo	Sorbic acid	12038012	1
Caelo	Sorbic acid	10120303	1
Caelo	Sorbic acid	12038025	1
Caelo	Sorbic acid	16208007	1
Caelo	Sorbic acid	16208015	9
Caelo	Sorbic acid	105505	1
Caelo	Sorbic acid	12038020	5
Caelo	Sorbic acid	12038002	1
Caelo	Sorbic acid	12038016	8
Caelo	Sorbic acid	16208018	2
Caelo	Sorbitan stearate	17178102	9
Caelo	Sorbitan stearate	17030805	4
Caelo	Sorbitan stearate	13B19-N01	2
Caelo	Sorbitan stearate	14218804	1
Caelo	Sorbitan stearate	15300302	7
Caelo	Sorbitan stearate	16235202	2
Caelo	Sorbitan stearate	16235207	9
Caelo	Sorbitan stearate	14038504	1
Caelo	Sorbitan stearate	14218808	3
Caelo	Sorbitan stearate	26031507	1
Caelo	Sorbitan stearate	12274813	2
Caelo	Sorbitan stearate	12274807	1
Caelo	Sorbitan stearate	14218803	3
Caelo	Sorbitan stearate	14038511	2
Caelo	Sorbitan stearate	14038510	1
Caelo	Sorbitan stearate	13209418	3
Caelo	Sorbitan stearate	13209409	3
Caelo	Sorbitan stearate	13209417	3
Caelo	Sorbitan stearate	18204705	1
Caelo	Sorbitan stearate	15415202	8
Caelo	Sorbitan stearate	13209419	3
Caelo	Sorbitan stearate	17345402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sorbitan stearate	15300308	4
Caelo	Sorbitan stearate	13209410	1
Caelo	Sorbitol	17366301	3
Caelo	Sorbitol	12249139	1
Caelo	Sorbitol	19017304	1
Caelo	Sorbitol	18163607	1
Caelo	Sorbitol	18163606	2
Caelo	Sorbitol	11379320	1
Caelo	Sorbitol	17366304	1
Caelo	Sorbitol	16101303	1
Caelo	Sorbitol	11379333	1
Caelo	Sorbitol	1714090903	1
Caelo	Sorbitol	17140903	6
Caelo	Sorbitol	113793	1
Caelo	Sorbitol	16231208	1
Caelo	Sorbitol	16353107	1
Caelo	Sorbitol	16038514	5
Caelo	Sorbitol	17140902	2
Caelo	Sorbitol	16353102	3
Caelo	Sorbitol	12249126	2
Caelo	Sorbitol	12249124	2
Caelo	Sorbitol	10383323	1
Caelo	Sorbitol	12345	1
Caelo	Sorbitol	23071201	1
Caelo	Sorbitol	17011305	1
Caelo	Sorbitol	11379332	2
Caelo	Sorbitol	11379337	2
Caelo	Sorbitol	12L07-N12	1
Caelo	Sorbitol	5021307	1
Caelo	Sorbitol	1083330	1
Caelo	Sorbitol	12249108	3
Caelo	Sorbitol	11379322	1
Caelo	Sorbitol	12249109	2
Caelo	Sorbitol	12249110	4
Caelo	Sorbitol	11379336	2
Caelo	Sorbitol	11379321	1
Caelo	Sorbitol	16038515	6
Caelo	Sorbitol	225282222169	1
Caelo	Sorbitol	15118912	1
Caelo	Sorbitol	15118902	2
Caelo	Sorbitol	14250806	1
Caelo	Sorbitol	26071402	1
Caelo	Sorbitol	8081402	1
Caelo	Sorbitol	23081401	1
Caelo	Sorbitol	4091401	1
Caelo	Sorbitol	15118901	3
Caelo	Sorbitol	14250808	1
Caelo	Sorbitol	14250807	1
Caelo	Sorbitol	20021501	1
Caelo	Sorbitol	13278241	1
Caelo	Sorbitol	13278218	1
Caelo	Sorbitol	13278247	1
Caelo	Sorbitol	13278239	1
Caelo	Sorbitol	13278250	1
Caelo	Sorbitol	13278216	2
Caelo	Sorbitol	133335307	1
Caelo	Sorbitol	15208302	4
Caelo	Sorbitol	16038503	1
Caelo	Sorbitol	15208310	4
Caelo	Sorbitol	16038502	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sorbitol	12249140	1
Caelo	Sorbitol	16038501	4
Caelo	Sorbitol	13335304	1
Caelo	Sorbitol	13335307	3
Caelo	Sorbitol	15208311	3
Caelo	Sorbitol	1512016-01	4
Caelo	Sorbitol	13278204	1
Caelo	Sorbitol	15208309	3
Caelo	Sorbitol	15208305	3
Caelo	Sorbitol	15012708	1
Caelo	Sorbitol	15208306	2
Caelo	Sorbitol	15208303	1
Caelo	Sorbitol		2
Caelo	Sorbitol	13278246	2
Caelo	Sorbitol	13278249	2
Caelo	Sorrel herb, cut	13289502	2
Caelo	Sorrel herb, cut	130546	1
Caelo	Sour cherry stems, cut	12042307	1
Caelo	Sour cherry stems, cut	10367019	1
Caelo	Sour cherry stems, cut	12042301	1
Caelo	Sour cherry stems, cut	10367014	1
Caelo	Soy lecithin (instant...	11371504	1
Caelo	Soy lecithin (instant...	14281618	4
Caelo	Soy lecithin (instant...	15262701	4
Caelo	Soy lecithin (instant...	12139011	3
Caelo	Soy lecithin (instant...	14196305	1
Caelo	Soy lecithin (instant...	13279814	1
Caelo	Soy lecithin (instant...	14281605	1
Caelo	Soy lecithin (instant...	14281606	1
Caelo	Soy lecithin (instant...	14281620	1
Caelo	Soy lecithin (instant...	13030503	2
Caelo	Soy lecithin (instant...	17081501	1
Caelo	Soy lecithin (instant...	12139001	3
Caelo	Soy lecithin (instant...	15262704	4
Caelo	Soy lecithin (instant...	17118101	1
Caelo	Soy lecithin (instant...	17118108	1
Caelo	Soy lecithin (instant...	17322401	1
Caelo	Soy lecithin (instant...	Nachprüfung	1
Caelo	Soy lecithin (instant...	12139012	1
Caelo	Soy lecithin (instant...	1526704	1
Caelo	Soy lecithin (instant...	13279806	1
Caelo	Spearmint leaf, cut	13105003	1
Caelo	Spearmint leaf, cut	12376706	1
Caelo	Spearmint leaf, cut	13105001	10
Caelo	Spearmint leaf, cut	12133901	3
Caelo	Spearmint leaf, cut	12133909	1
Caelo	Spearmint leaf, cut	1512Q-01505	1
Caelo	Spearmint leaf, cut	12133906	2
Caelo	Spearmint leaf, cut	11037409+110374	1
Caelo	Spironolactone	13124304	1
Caelo	Spironolactone	14378717	2
Caelo	Spironolactone	13190309	1
Caelo	Spironolactone	13124314	1
Caelo	Spironolactone	1314309	1
Caelo	Spironolactone	13124309	3
Caelo	Spironolactone	1437807	1
Caelo	Spironolactone	15G22-B10	1
Caelo	Spironolactone	14378715	5
Caelo	Spironolactone	14378721	3
Caelo	Spironolactone	14378718	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Spironolactone	13124304	6
Caelo	Spironolactone	14378701	2
Caelo	Spironolactone	12049109	1
Caelo	Spironolactone	13124303	8
Caelo	Spironolactone	12049118	2
Caelo	Spironolactone	2013U-03249	1
Caelo	Spironolactone	11042602	1
Caelo	Spironolactone	11G28-N03	1
Caelo	Spironolactone		1
Caelo	Spironolactone	14378710	4
Caelo	Spironolactone	12049123	3
Caelo	Spironolactone	10159914	1
Caelo	Spironolactone	12049110	2
Caelo	Spironolactone	13124315	1
Caelo	Spironolactone	14378709	3
Caelo	Spironolactone	13124314	1
Caelo	Spironolactone	13124321	1
Caelo	Spironolactone	14378703	1
Caelo	Spironolactone	11042607	1
Caelo	Spironolactone	18208010	2
Caelo	Spironolactone	18208008	2
Caelo	Spironolactone	16077302	10
Caelo	Spironolactone	13124321	1
Caelo	Spironolactone	16077319	4
Caelo	Spironolactone	17116501	1
Caelo	Spironolactone	16A14-B05-317123	1
Caelo	Spironolactone	13124317	3
Caelo	Spironolactone	16077307	1
Caelo	Spironolactone	16077320	6
Caelo	Spironolactone	16077312	11
Caelo	Spironolactone	1438721	1
Caelo	Spironolactone	5013Q-03249	1
Caelo	Spironolactone	16H20-B06-322016	1
Caelo	Spironolactone	14378719	1
Caelo	Spironolactone	18208003	2
Caelo	Spironolactone	16077311	10
Caelo	Spironolactone	13124320	3
Caelo	Spironolactone	16077309	5
Caelo	Spironolactone	17116502	7
Caelo	Spironolactone	171165012	1
Caelo	Spironolactone	18004920	5
Caelo	Spironolactone	16077303	12
Caelo	Spironolactone	160773	1
Caelo	Spironolactone	16077306	3
Caelo	Spironolactone	10H24-N10	1
Caelo	Spironolactone	18004915	3
Caelo	Spironolactone	140202901	1
Caelo	Spironolactone	18004917	10
Caelo	Spironolactone	18004904	8
Caelo	Spironolactone	18004803	7
Caelo	Spironolactone	16077325	8
Caelo	Spironolactone	16077305	5
Caelo	Spironolactone	14378702	1
Caelo	Spironolactone	16077326	6
Caelo	Spironolactone	16077324	1
Caelo	Spironolactone	18004908	3
Caelo	Spray-dried arabic gum	3314A160121	1
Caelo	Spray-dried arabic gum	11345715	1
Caelo	Spray-dried arabic gum	12138607	1
Caelo	Spray-dried arabic gum	15413207	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Spray-dried arabic gum	12372302	1
Caelo	Spray-dried arabic gum	12138612	1
Caelo	Spray-dried arabic gum	12111301	1
Caelo	Spray-dried arabic gum	15086704	1
Caelo	Spray-dried arabic gum	15413201	1
Caelo	St. John's wort, cut	111185608	1
Caelo	St. John's wort, cut	12092209	2
Caelo	St. John's wort, cut	13091303	1
Caelo	St. John's wort, cut	12092204	1
Caelo	St. John's wort, cut	16041301	1
Caelo	St. John's wort, cut	11185607	1
Caelo	St. John's wort, cut	13203204	6
Caelo	St. John's wort, cut	1185608	1
Caelo	St. John's wort, cut	12092201	1
Caelo	St. John's wort, cut	4021B110321	1
Caelo	St. John's wort, cut	15068105	1
Caelo	St. John's wort, cut	132203204	1
Caelo	St. John's wort, cut	11185602	1
Caelo	St. John's wort, cut	11185608	18
Caelo	St. John's wort, cut	30011401	1
Caelo	St. John's wort, cut	5021301	1
Caelo	St. John's wort, cut	13203205	1
Caelo	St. John's wort, cut	268913	1
Caelo	Star anise	10169101	1
Caelo	Star anise	10376405	1
Caelo	Star anise	11314104	1
Caelo	Star anise	12226211	1
Caelo	Star anise	11314103	1
Caelo	Star anise	10376403	2
Caelo	Star anise, cut	16253201	2
Caelo	Stearic acid 50	18080901	1
Caelo	Stearic acid 50	17081006	3
Caelo	Stearic acid 50	17081003	1
Caelo	Stearic acid 50	1617081006	1
Caelo	Stearic acid 50	14359314	3
Caelo	Stearic acid 50	11308505	2
Caelo	Stearic acid 50	11103706	1
Caelo	Stearic acid 50	1609611	1
Caelo	Stearic acid 50	11308514	1
Caelo	Stearic acid 50	12241911	2
Caelo	Stearic acid 50	12241915	3
Caelo	Stearic acid 50	12241914	1
Caelo	Stearic acid 50	12241930	3
Caelo	Stearic acid 50	14359314	1
Caelo	Stearic acid 50	16096311	1
Caelo	Stearic acid 50	16096301	1
Caelo	Stemless carline root	25111306	1
Caelo	Stemless carline root	13033501	1
Caelo	Stemless carline root	13130304	1
Caelo	Stemless carline root	14347401	1
Caelo	Stemless carline root	13130303	1
Caelo	Stemless carline root	13130301	6
Caelo	Stemless carline root	23071203	1
Caelo	Stemless carline root	12232803	1
Caelo	Stemless carline root	12144804	1
Caelo	Strawberry leaf, cut	11117909	1
Caelo	Strawberry leaf, cut	5122013c	1
Caelo	Sucrose	14340006	1
Caelo	Sucrose	15182804	23
Caelo	Sucrose	162061	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sucrose	15182805	4
Caelo	Sucrose	15446007	1
Caelo	Sucrose	16328605	17
Caelo	Sucrose	16051305	1
Caelo	Sucrose	16328606	3
Caelo	Sucrose	16206103	5
Caelo	Sucrose	15446002	11
Caelo	Sucrose	15446003	6
Caelo	Sucrose	13352916	1
Caelo	Sucrose	16049902	2
Caelo	Sucrose	16051301	1
Caelo	Sucrose	17354506	6
Caelo	Sucrose	16328603	8
Caelo	Sucrose	13352902	2
Caelo	Sucrose	13352916	2
Caelo	Sucrose	16051304	1
Caelo	Sucrose	160499	1
Caelo	Sucrose	17021620306	1
Caelo	Sucrose	14197302	2
Caelo	Sucrose	1352920	1
Caelo	Sucrose	14197313	1
Caelo	Sucrose	28041504	1
Caelo	Sucrose	15028801	11
Caelo	Sucrose	14340004	3
Caelo	Sucrose	14340008	1
Caelo	Sucrose	14340009	4
Caelo	Sucrose	20031515	1
Caelo	Sucrose	14197301	3
Caelo	Sucrose	1802445	1
Caelo	Sucrose	18154603	11
Caelo	Sucrose	13352920	1
Caelo	Sucrose	14340003	2
Caelo	Sucrose	1430009	1
Caelo	Sucrose	14340003	1
Caelo	Sucrose	27111802	1
Caelo	Sucrose	24061502	1
Caelo	Sucrose	18257402	16
Caelo	Sucrose	18281803	1
Caelo	Sucrose	18257401	1
Caelo	Sucrose	17288702	1
Caelo	Sucrose	18281804	1
Caelo	Sucrose	18329702	1
Caelo	Sucrose	14197301	1
Caelo	Sucrose	29081502	1
Caelo	Sucrose	15028802	5
Caelo	Sucrose	8024405	1
Caelo	Sucrose	14197302	1
Caelo	Sucrose	13352920	6
Caelo	Sucrose	12107806	2
Caelo	Sucrose	12241804	2
Caelo	Sucrose	17106601	1
Caelo	Sucrose	12241808	5
Caelo	Sucrose	12241803	2
Caelo	Sucrose	13031417	2
Caelo	Sucrose	130031417	1
Caelo	Sucrose	12107805	1
Caelo	Sucrose	13063907	1
Caelo	Sucrose	13031410	2
Caelo	Sucrose	13031406	2
Caelo	Sucrose	13031414	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Sucrose	20011401	1
Caelo	Sucrose	17106604	8
Caelo	Sucrose	17106603	32
Caelo	Sucrose	13352911	1
Caelo	Sucrose	17042103	2
Caelo	Sucrose	13111706	1
Caelo	Sucrose	11246306	1
Caelo	Sucrose	1649902	1
Caelo	Sucrose	12241801	1
Caelo	Sucrose	13352902	10
Caelo	Sucrose	12241817	2
Caelo	Sucrose	12280309	1
Caelo	Sucrose	13352909	4
Caelo	Sucrose	13352905	1
Caelo	Sucrose		1
Caelo	Sucrose	81172478	2
Caelo	Sucrose	16328607	11
Caelo	Sucrose	16049907	13
Caelo	Sucrose	17247203	1
Caelo	Sucrose	18024407	2
Caelo	Sucrose	18024405	11
Caelo	Sucrose	13031409	1
Caelo	Sucrose	13352908	2
Caelo	Sucrose	11246313	2
Caelo	Sucrose	13031405	3
Caelo	Sucrose	28081308	1
Caelo	Sucrose	13352901	2
Caelo	Sucrose	154460	1
Caelo	Sucrose	15156008	1
Caelo	Sulfacetamide	11095912	1
Caelo	Sulfacetamide	11085912	4
Caelo	Sulfacetamide	10145811	1
Caelo	Sulfacetamide	11085911	4
Caelo	Sulfacetamide	11085909	1
Caelo	Sulfacetamide	11085913	1
Caelo	Sulfacetamide	15445901	1
Caelo	Sulfacetamide	10145802	1
Caelo	Sulfacetamide	11085914	2
Caelo	Sulfadiazine	17358910	3
Caelo	Sulfadiazine	17358905	2
Caelo	Sulfadiazine	17358901	1
Caelo	Sulfadiazine	13150004	1
Caelo	Sulfadiazine	11194204	1
Caelo	Sulfamerazine	10218401	1
Caelo	Sulfogaiacol hydrate	15343912	12
Caelo	Sweet woodruff herb, ...	12103509	1
Caelo	Sweet woodruff herb, ...	12103513	2
Caelo	Sweet woodruff herb, ...	12103514	1
Caelo	Sweet woodruff herb, ...	10375802	1
Caelo	Sweet woodruff herb, ...	52236077	1
Caelo	Sweet woodruff herb, ...	10375812	1
Caelo	Sweet woodruff herb, ...	4409M-01621	1
Caelo	Sweet woodruff herb, ...	18031304	1
Caelo	Sweet woodruff herb, ...	12103508	1
Caelo	Sweet woodruff herb, ...	12103510	1
Caelo	SyrSpend [®] SF Alka	18871B06	1
Caelo	SyrSpend [®] SF Alka	18J25-B15-357803	1
Caelo	SyrSpend [®] SF Alka	353473	1
Caelo	SyrSpend [®] SF Alka	18C13-B11-350851	2
Caelo	SyrSpend [®] SF Alka	720026-000	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Talcum	12267401	1
Caelo	Talcum	13146308	1
Caelo	Talcum	13146307	1
Caelo	Tall golden rod, cut	13386901	1
Caelo	Tall golden rod, cut	13386903	1
Caelo	Tall golden rod, cut	14018404	1
Caelo	Tall golden rod, cut	12343109	1
Caelo	Tall golden rod, cut	4021405	1
Caelo	Tall golden rod, cut	28051302	1
Caelo	Tall golden rod, cut	12343105	1
Caelo	Tall golden rod, cut	27121302	1
Caelo	Tall golden rod, cut	6081302	1
Caelo	Tall golden rod, cut	11207010	3
Caelo	Tall golden rod, cut	15061203	1
Caelo	Tall golden rod, cut	265283	1
Caelo	Tall golden rod, cut	11280301	4
Caelo	Tall golden rod, cut	17091203	1
Caelo	Tall golden rod, cut	11280308	9
Caelo	Tall golden rod, cut	20111209	1
Caelo	Tall golden rod, cut	112080308	1
Caelo	Tall golden rod, cut	12336801	1
Caelo	Tall golden rod, cut	11280310	6
Caelo	Tall golden rod, cut	11280309	1
Caelo	Tall golden rod, cut	12319003	22
Caelo	Tall golden rod, cut	12319002	2
Caelo	Tall golden rod, cut	12343106	1
Caelo	Tall golden rod, cut	14112012	1
Caelo	Tall golden rod, cut	12343103	17
Caelo	Tall golden rod, cut	2656153	1
Caelo	Tall golden rod, cut	4113A-01684	1
Caelo	Tall golden rod, cut	12319004	10
Caelo	Tall golden rod, cut	12319006	3
Caelo	Tall golden rod, cut	12343101	3
Caelo	Tall golden rod, cut	12319007	21
Caelo	Tall golden rod, cut	12343102	3
Caelo	Tall golden rod, cut	314050-1-p	1
Caelo	Tall golden rod, cut	4091302	1
Caelo	Tall golden rod, cut	12343104	5
Caelo	Tall hawkweed, cut	3611b111005	1
Caelo	Tall nasturtium herb,...	13104401	4
Caelo	Tall nasturtium herb,...	107063	1
Caelo	Tall nasturtium herb,...	12031602	1
Caelo	Tall nasturtium herb,...	13104402	2
Caelo	Tall nasturtium herb,...	12031605	1
Caelo	Tall nasturtium herb,...	12031601	1
Caelo	Tall nasturtium herb,...	12031604	2
Caelo	Tall nasturtium herb,...	13104403	5
Caelo	Tall nasturtium herb,...	12031603	1
Caelo	Tall nasturtium herb,...	13104404	2
Caelo	Tall nasturtium herb,...	13104405	1
Caelo	Tall nasturtium herb,...	13104406	1
Caelo	Tall nasturtium herb,...	12031606	1
Caelo	Tannin	1524308	1
Caelo	Tannin	15091506	1
Caelo	Tannin	15243802	3
Caelo	Tannin	15091503	2
Caelo	Tannin	10019007	1
Caelo	Tannin	15091504	4
Caelo	Tannin	14031913	2
Caelo	Tannin	14031918	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tannin	14031915	2
Caelo	Tannin	28011502	1
Caelo	Tannin	14031912	5
Caelo	Tannin	14031905	3
Caelo	Tannin	14031903	5
Caelo	Tannin	14031901	2
Caelo	Tannin	13031401	1
Caelo	Tannin	12336014	6
Caelo	Tannin	12336013	6
Caelo	Tannin	2211Q-01052	1
Caelo	Tannin	15243801	9
Caelo	Tannin	12336005	2
Caelo	Tannin	1236001	1
Caelo	Tannin	12336009	10
Caelo	Tannin	21101305	1
Caelo	Tannin	12336001	12
Caelo	Tannin	12171106	8
Caelo	Tannin	12171103	2
Caelo	Tannin	9128608	1
Caelo	Tannin	10315110	3
Caelo	Tannin	12171101	4
Caelo	Tannin	10315107	2
Caelo	Tannin	14031919	6
Caelo	Tannin	14031920	1
Caelo	Tannin	12336007	2
Caelo	Tannin	11302307	6
Caelo	Tannin	4413I1052	1
Caelo	Tannin	15091509	17
Caelo	Tannin	17343515	1
Caelo	Tannin	17343505	1
Caelo	Tannin	17343509	4
Caelo	Tannin	1509509	1
Caelo	Tannin	17343504	2
Caelo	Tannin	18017701	3
Caelo	Tannin	17343503	1
Caelo	Tannin		1
Caelo	Tannin	16074518	6
Caelo	Tannin	17343502	1
Caelo	Tannin	16074516	5
Caelo	Tannin	1217112263904	1
Caelo	Tannin	16074509/171027T	1
Caelo	Tannin	16D08-B04-328410	2
Caelo	Tannin	16074505	1
Caelo	Tannin	16074504	13
Caelo	Tannin	15091515	10
Caelo	Tannin	16074517	1
Caelo	Tannin	15091511	2
Caelo	Tannin	15091512	1
Caelo	Tartaric acid	15349106	2
Caelo	Tartaric acid	13125-B01-294093	1
Caelo	Tartaric acid	10180713	2
Caelo	Tartaric acid	1169406	1
Caelo	Tartaric acid	15164010	1
Caelo	Tartaric acid	15164008	3
Caelo	Tartaric acid	12044103	5
Caelo	Tartaric acid	15130206	1
Caelo	Tartaric acid	15349109	4
Caelo	Tartaric acid	15349101	3
Caelo	Tartaric acid	12240413	1
Caelo	Tartaric acid	14G07-B02-308555	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tartaric acid	10199713	1
Caelo	Tartaric acid	15349105	6
Caelo	Tartaric acid	151302	1
Caelo	Tartaric acid	12286904	3
Caelo	Tartaric acid	15130207	1
Caelo	Tartaric acid	15164002	3
Caelo	Tartaric acid	12286932	1
Caelo	Tartaric acid	15130201	6
Caelo	Tartaric acid	12286911	1
Caelo	Tartaric acid	10180708	2
Caelo	Tartaric acid	13116306	1
Caelo	Tartaric acid	12044101	2
Caelo	Tartaric acid	11k10-n07	1
Caelo	Tartaric acid	11169406	2
Caelo	Tartaric acid	12286924	1
Caelo	Tartaric acid	17091402	1
Caelo	Tartaric acid	15130205	2
Caelo	Tartaric acid	13217807	1
Caelo	Tartaric acid	9041703	1
Caelo	Tartaric acid	14054521	1
Caelo	Tartaric acid	13217806	2
Caelo	Tartaric acid	13116307	1
Caelo	Tartaric acid	12286939	1
Caelo	Tartaric acid	14186213	4
Caelo	Tartaric acid	15045304	1
Caelo	Tartaric acid	15164001	4
Caelo	Tartaric acid	12286934	1
Caelo	Tartaric acid	15349102	2
Caelo	Tartaric acid	16260404	5
Caelo	Tartaric acid	16260403	2
Caelo	Tartaric acid	17337504	1
Caelo	Tartaric acid	18121208	1
Caelo	Tartaric acid	18121210	1
Caelo	Tartaric acid	18121201	1
Caelo	Tartaric acid	1812104	1
Caelo	Tartaric acid	9494239	1
Caelo	Tartaric acid	17337505	1
Caelo	Tartaric acid	12286906	5
Caelo	Tartaric acid	1204102	1
Caelo	Tartaric acid	12H02-n07	1
Caelo	Tartaric acid	4110N-01053	1
Caelo	Tartaric acid	12286903	2
Caelo	Tartaric acid	10180714	1
Caelo	Tartaric acid	12044119	3
Caelo	Tartaric acid	12044102	3
Caelo	Tartaric acid	12240402	1
Caelo	Tartaric acid	18121211	1
Caelo	Tartaric acid	15349110	2
Caelo	Tartaric acid	16260406	4
Caelo	Tartaric acid	17347103	3
Caelo	Tartaric acid	17347104	3
Caelo	Tartaric acid	16260405	1
Caelo	Tartaric acid	12240415	1
Caelo	Tartaric acid	8051809	1
Caelo	Testosterone	12176305	3
Caelo	Testosterone	15412202	2
Caelo	Testosterone	13406004	2
Caelo	Testosterone	13134804	1
Caelo	Testosterone	18147803	2
Caelo	Testosterone	APL00034466T	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Testosterone	14355605	1
Caelo	Testosterone	AP5106061T	6
Caelo	Testosterone	17119911	1
Caelo	Testosterone	17119902	1
Caelo	Testosterone	14201105	1
Caelo	Testosterone	13134806	4
Caelo	Testosterone	15195803	1
Caelo	Testosterone	12176301	1
Caelo	Testosterone	15195802	1
Caelo	Testosterone	12176311Mg: 1g	1
Caelo	Testosterone	15412201	3
Caelo	Testosterone	11/07-N12	1
Caelo	Testosterone propionate	1263607	1
Caelo	Testosterone propionate	36814m1	1
Caelo	Testosterone propionate	12163617	5
Caelo	Testosterone propionate	12163619	2
Caelo	Testosterone propionate	36814M1	3
Caelo	Testosterone propionate	15169504	1
Caelo	Testosterone propionate	15003002	3
Caelo	Testosterone propionate	12163616	6
Caelo	Testosterone propionate	15444703	12
Caelo	Testosterone propionate	15003004	3
Caelo	Testosterone propionate	614E2376	1
Caelo	Testosterone propionate	12163613	2
Caelo	Testosterone propionate	15003008	2
Caelo	Testosterone propionate	15003002	5
Caelo	Testosterone propionate	103004	1
Caelo	Testosterone propionate	12163621	1
Caelo	Testosterone propionate	15444705	3
Caelo	Testosterone propionate	17325707	5
Caelo	Testosterone propionate	15444707	4
Caelo	Testosterone propionate	106212376	1
Caelo	Testosterone propionate	1108212376	1
Caelo	Testosterone propionate	12060102	2
Caelo	Testosterone propionate	612212376	1
Caelo	Testosterone propionate	151100301	1
Caelo	Testosterone propionate	903312376	1
Caelo	Testosterone propionate	10235809	1
Caelo	Testosterone propionate	12163615	2
Caelo	Testosterone propionate	17325709	1
Caelo	Testosterone propionate	17325719	1
Caelo	Testosterone propionate	2013U-02376	1
Caelo	Testosterone propionate	18B14-B02-353456	1
Caelo	Testosterone propionate	17325712	2
Caelo	Testosterone propionate	18b06-B04	1
Caelo	Testosterone propionate	12163605	2
Caelo	Testosterone propionate	15444712	1
Caelo	Testosterone propionate	15444721	1
Caelo	Testosterone propionate	15444709	1
Caelo	Testosterone propionate	17J16-B03-341320	1
Caelo	Testosterone propionate	15444723	1
Caelo	Testosterone propionate	1710027-01	3
Caelo	Testosterone propionate	15444717	1
Caelo	Testosterone propionate	15444714	2
Caelo	Testosterone propionate	15444719	1
Caelo	Testosterone propionate	15444710	5
Caelo	Testosterone propionate	15444711	1
Caelo	Testosterone propionate	12163612	6
Caelo	Testosterone propionate	15444706	4
Caelo	Tetracaine base	T-290413-1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tetracaine base	16E23-B04	1
Caelo	Tetracaine base	17111706	5
Caelo	Tetracaine base	1516B06312783	1
Caelo	Tetracaine base	16A13-B07-318736	9
Caelo	Tetracaine base	14A16-B04-292916	3
Caelo	Tetracaine base	16a13-b07-318736	1
Caelo	Tetracaine base	14A16-B04	1
Caelo	Tetracaine base	14A16-B04-290207	1
Caelo	Tetracaine base	17111705	1
Caelo	Tetracaine base	1805019-01	2
Caelo	Tetracaine hydrochlor...	12019802	4
Caelo	Tetracaine hydrochlor...	10118211	1
Caelo	Tetracaine hydrochlor...	10118210	2
Caelo	Tetracaine hydrochlor...	12019812	3
Caelo	Tetracaine hydrochlor...	12019810	5
Caelo	Tetracaine hydrochlor...	18060609	1
Caelo	Tetracaine hydrochlor...	18060602	4
Caelo	Tetracaine hydrochlor...	12019816	12
Caelo	Tetracaine hydrochlor...	18060621	1
Caelo	Tetracaine hydrochlor...	18060612	3
Caelo	Tetracaine hydrochlor...	12019815	1
Caelo	Tetracaine hydrochlor...	12019822	6
Caelo	Tetracaine hydrochlor...	9070003	2
Caelo	Tetracaine hydrochlor...	10118225	4
Caelo	Tetracaine hydrochlor...	12019820	1
Caelo	Tetracaine hydrochlor...	16G28-B03-332099	2
Caelo	Tetracaine hydrochlor...	16G25-B03-332099	1
Caelo	Tetracaine hydrochlor...	15437119	6
Caelo	Tetracaine hydrochlor...	15437114	5
Caelo	Tetracaine hydrochlor...	1110010-07	1
Caelo	Tetracaine hydrochlor...	15K04-B03-318369	2
Caelo	Tetracaine hydrochlor...	1537109	1
Caelo	Tetracaine hydrochlor...	15437108	1
Caelo	Tetracaine hydrochlor...	13155705	3
Caelo	Tetracaine hydrochlor...	15437109	9
Caelo	Tetracaine hydrochlor...	13155711	10
Caelo	Tetracaine hydrochlor...	150604002	1
Caelo	Tetracaine hydrochlor...	15437106	2
Caelo	Tetracaine hydrochlor...	15437102	1
Caelo	Tetracaine hydrochlor...	15356902	2
Caelo	Tetracaine hydrochlor...	15196104	2
Caelo	Tetracaine hydrochlor...	15079507	1
Caelo	Tetracaine hydrochlor...	15H24-B03-311450	7
Caelo	Tetracaine hydrochlor...	6093549	3
Caelo	Tetracaine hydrochlor...	14G14-B08	2
Caelo	Tetracaine hydrochlor...	170601801	1
Caelo	Tetracaine hydrochlor...	15196103	5
Caelo	Tetracaine hydrochlor...	13155724	4
Caelo	Tetracaine hydrochlor...	15196102	4
Caelo	Tetracaine hydrochlor...	13155714	3
Caelo	Tetracaine hydrochlor...	14C18-B03-300200	5
Caelo	Tetracaine hydrochlor...	13155727	6
Caelo	Tetracaine hydrochlor...	13155703	1
Caelo	Tetracaine hydrochlor...	1508040-01	1
Caelo	Tetracaine hydrochlor...	1260917/15437119	1
Caelo	Tetracaine hydrochlor...	259/22	1
Caelo	Tetracaine hydrochlor...	1313012-01	1
Caelo	Tetracaine hydrochlor...	13155706	11
Caelo	Tetracaine hydrochlor...	14c18-b03-293966	1
Caelo	Tetracaine hydrochlor...	15437122	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tetracaine hydrochlor...	13155716	2
Caelo	Tetracaine hydrochlor...	15437121	4
Caelo	Tetracaine hydrochlor...	15079501	2
Caelo	Tetracaine hydrochlor...	27111305	1
Caelo	Tetracaine hydrochlor...	13155720	2
Caelo	Tetracycline hydrochl...	1017A-02378	1
Caelo	Tetracycline hydrochl...	16116706	5
Caelo	Tetracycline hydrochl...	15230711	1
Caelo	Tetracycline hydrochl...	12353602	1
Caelo	Tetracycline hydrochl...	1160117	1
Caelo	Tetracycline hydrochl...	16116701	3
Caelo	Tetracycline hydrochl...	14364601	1
Caelo	Tetracycline hydrochl...	16116704	1
Caelo	Tetracycline hydrochl...	15230715	2
Caelo	Tetracycline hydrochl...	15230720	2
Caelo	Tetracycline hydrochl...	15230717	3
Caelo	Tetracycline hydrochl...	4315I-02378	1
Caelo	Tetracycline hydrochl...	15230707	2
Caelo	Tetracycline hydrochl...	15230714	1
Caelo	Tetracycline hydrochl...	14080605	2
Caelo	Tetracycline hydrochl...	2714U-02378	1
Caelo	Tetracycline hydrochl...	16116720	1
Caelo	Tetracycline hydrochl...	12353612	1
Caelo	Tetracycline hydrochl...	14364612	2
Caelo	Tetracycline hydrochl...	14364614	1
Caelo	Tetracycline hydrochl...	14364616	2
Caelo	Tetracycline hydrochl...	14H06-B05-313591	1
Caelo	Tetracycline hydrochl...	14364604	4
Caelo	Tetracycline hydrochl...	12353609	1
Caelo	Tetracycline hydrochl...	12353603	6
Caelo	Tetracycline hydrochl...	12138104	1
Caelo	Tetracycline hydrochl...	166116709	1
Caelo	Tetracycline hydrochl...	18308302	1
Caelo	Tetracycline hydrochl...	4012Q-02378	1
Caelo	Tetracycline hydrochl...	18308303	1
Caelo	Tetracycline hydrochl...	103201301	1
Caelo	Tetracycline hydrochl...	13201302	1
Caelo	Tetracycline hydrochl...	1402030-01	10
Caelo	Tetracycline hydrochl...	12353613	1
Caelo	Tetracycline hydrochl...	18198608	2
Caelo	Tetracycline hydrochl...	13L10B02296956	1
Caelo	Tetracycline hydrochl...	18198602	4
Caelo	Tetracycline hydrochl...	17285808	3
Caelo	Tetracycline hydrochl...	17285809	4
Caelo	Tetracycline hydrochl...	172855809	1
Caelo	Tetracycline hydrochl...	1611617	1
Caelo	Tetracycline hydrochl...	15230703	4
Caelo	Tetracycline hydrochl...	17285803	4
Caelo	Tetracycline hydrochl...	5091202	1
Caelo	Tetracycline hydrochl...	111002901	1
Caelo	Tetracycline hydrochl...	9248704	1
Caelo	Tetracycline hydrochl...	16116721	3
Caelo	Tetracycline hydrochl...	13201301	4
Caelo	Tetracycline hydrochl...	1204039	2
Caelo	Tetracycline hydrochl...	17285801	4
Caelo	Tetracycline hydrochl...	12138106	2
Caelo	Tetracycline hydrochl...	13201309	1
Caelo	Tetracycline hydrochl...	16116714	1
Caelo	Tetracycline hydrochl...	16116715	5
Caelo	Tetracycline hydrochl...	16116709	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tetracycline hydrochl. . .	14080606	3
Caelo	Theophylline	12166305	3
Caelo	Theophylline	17196203	1
Caelo	Theophylline	15203010	1
Caelo	Theophylline	12341104	1
Caelo	Theophylline	12341101	2
Caelo	Theophylline	12166306	1
Caelo	Theophylline	12166311	1
Caelo	Theophylline	12166312	3
Caelo	Theophylline-ethylene. . .	15380302	1
Caelo	Theriac, powdered	14156501	1
Caelo	Theriac, powdered	12289305	2
Caelo	Theriac, powdered	11082413	1
Caelo	Theriac, powdered	15312508	1
Caelo	Theriac, powdered	15312511	1
Caelo	Theriac, powdered	12289307	8
Caelo	Theriac, powdered	15101307	1
Caelo	Theriac, powdered	12189404	2
Caelo	Theriac, powdered	15312504	1
Caelo	Theriac, powdered	412091-2P	1
Caelo	Theriac, powdered	12289303	1
Caelo	Theriac, powdered	1531501	1
Caelo	Theriac, powdered	2021301	1
Caelo	Theriac, powdered	12289312	1
Caelo	Theriac, powdered	15312501	1
Caelo	Theriac, powdered	11081515	1
Caelo	Theriac, powdered	15026401	1
Caelo	Theriac, powdered	13407304	2
Caelo	Theriac, powdered	12289301	4
Caelo	Theriac, powdered	12189401	1
Caelo	Theriac, powdered	12289306	1
Caelo	Thiamine chloride hyd. . .	11371319	3
Caelo	Thiamine chloride hyd. . .		1
Caelo	Thiamine chloride hyd. . .	11371309	2
Caelo	Thiamine chloride hyd. . .	11061401	1
Caelo	Thiamine chloride hyd. . .	4071402	1
Caelo	Thiamine chloride hyd. . .	12340712	2
Caelo	Thiamine chloride hyd. . .	14069502	1
Caelo	Thiamine chloride hyd. . .	12091407	1
Caelo	Thiamine chloride hyd. . .	11371313	1
Caelo	Thiamine chloride hyd. . .	8121405	1
Caelo	Thiamine chloride hyd. . .	13296610	4
Caelo	Thiamine chloride hyd. . .	15003201	1
Caelo	Thiamine chloride hyd. . .	15199009	2
Caelo	Thiamine chloride hyd. . .	5122	1
Caelo	Thiamine chloride hyd. . .	16015609	1
Caelo	Thiamine chloride hyd. . .	15199099	1
Caelo	Thiamine chloride hyd. . .	171201301	1
Caelo	Thiamine chloride hyd. . .	15269613	1
Caelo	Thiamine chloride hyd. . .	16015607	1
Caelo	Thiomersal	10296413	1
Caelo	Thyme, rubbed	13011402	1
Caelo	Thyme, rubbed	15101201	1
Caelo	Thyme, rubbed		1
Caelo	Thyme, rubbed	16011401	1
Caelo	Thyme, rubbed	11345007	5
Caelo	Thyme, rubbed	11121306	1
Caelo	Thyme, rubbed	13048309-3475312161	1
Caelo	Thyme, rubbed	11345011	12
Caelo	Thyme, rubbed	13048309	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Thyme, rubbed	4021512	1
Caelo	Thyme, rubbed	26021301	1
Caelo	Thyme, rubbed	13048314	1
Caelo	Thyme, rubbed	22021301	1
Caelo	Thyme, rubbed	11021301	1
Caelo	Thyme, rubbed	29011303	1
Caelo	Thyme, rubbed	8061206	1
Caelo	Thyme, rubbed	12111201	1
Caelo	Thyme, rubbed	11345004	4
Caelo	Thyme, rubbed	14119001	1
Caelo	Thyme, rubbed	11345010	6
Caelo	Thyme, rubbed	11345001	1
Caelo	Thyme, rubbed	21121204	1
Caelo	Thyme, rubbed	15104401	1
Caelo	Thyme, rubbed	11345006	3
Caelo	Thyme, rubbed	23021304	1
Caelo	Thyme, rubbed	4041301	1
Caelo	Thyme, rubbed	25041302	1
Caelo	Thyme, rubbed	13048308	2
Caelo	Thyme, rubbed	25091304	1
Caelo	Thyme, rubbed	10365807	2
Caelo	Thyme, rubbed	13048301	14
Caelo	Thyme, rubbed	10365816	1
Caelo	Thyme, rubbed	13048302	2
Caelo	Thyme, rubbed	13081301	1
Caelo	Thyme, rubbed	29101301	1
Caelo	Thyme, rubbed	11345014	12
Caelo	Thyme, rubbed	11345002	1
Caelo	Thyme, rubbed	1134514	1
Caelo	Thyme, rubbed	11345013	9
Caelo	Thyme, rubbed	10365815	2
Caelo	Thymol	14025816	1
Caelo	Thymol	1730711	1
Caelo	Thymol	14025807	1
Caelo	Thymol	14025809	2
Caelo	Thymol	14302401	1
Caelo	Thymol	17324203	2
Caelo	Thymol	14310605	1
Caelo	Thymol	15130106	2
Caelo	Thymol	11178718	1
Caelo	Thymol	14025805	1
Caelo	Thymol	13147502	2
Caelo	Thymol	14025815	1
Caelo	Thymol	16144101	5
Caelo	Thymol	14310601	1
Caelo	Thymol	16144105	7
Caelo	Thymol	19061303	1
Caelo	Thymol	13288802	4
Caelo	Thymol	161441108	1
Caelo	Thymol	21042	1
Caelo	Thymol	11358616	1
Caelo	Thymol	2031503	1
Caelo	Thymol	15130101	1
Caelo	Thymol	11358608	1
Caelo	Thymol	11358607	2
Caelo	Thymol	15362908	3
Caelo	Thymol	15362901	6
Caelo	Thymol	16144108	3
Caelo	Tonka beans	20111312	1
Caelo	Tonka beans	12144101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tonka beans	12145101	2
Caelo	Tonka beans	1078C	1
Caelo	Tormentil	15085501	1
Caelo	Tormentil	15364302	4
Caelo	Tormentil	15178202	2
Caelo	Tormentil	15178203	1
Caelo	Tormentil	4753268	1
Caelo	Tormentil	14360103	1
Caelo	Tormentil	15178206	1
Caelo	Tormentil	1914A-02188	1
Caelo	Tormentil	12216803	4
Caelo	Tormentil	14008001	2
Caelo	Tormentil	14008002	2
Caelo	Tormentil	25022014A	1
Caelo	Tormentil	14008005	1
Caelo	Tormentil	11102010	1
Caelo	Tormentil	12040704	2
Caelo	Tormentil	10081203	1
Caelo	Tormentil	12040708	9
Caelo	Tormentil	2512A-02188	1
Caelo	Tormentil	14140001	1
Caelo	Tormentil	12040717	1
Caelo	Tormentil	12216805	1
Caelo	Tormentil	12216801	4
Caelo	Tormentil	12216802	1
Caelo	Tormentil	30011405	1
Caelo	Tormentil	12040713	1
Caelo	Tormentil	411082-2-P	1
Caelo	Tormentil	13220502	1
Caelo	Tormentil	12040707	1
Caelo	Tosylchloramide sodium	12277005	3
Caelo	Tosylchloramide sodium	15370501	4
Caelo	Tosylchloramide sodium	15370505	2
Caelo	Tosylchloramide sodium	16236603	5
Caelo	Tosylchloramide sodium	17061503	2
Caelo	Tosylchloramide sodium	16236603/1180117	1
Caelo	Tosylchloramide sodium	15183905	1
Caelo	Tosylchloramide sodium	17061511	1
Caelo	Tosylchloramide sodium	16236611	1
Caelo	Tosylchloramide sodium	17226503	2
Caelo	Tosylchloramide sodium	16236607	9
Caelo	Tosylchloramide sodium	15039603	4
Caelo	Tosylchloramide sodium	1310206	1
Caelo	Tosylchloramide sodium	13361607	6
Caelo	Tosylchloramide sodium	13361606	3
Caelo	Tosylchloramide sodium	13103210	2
Caelo	Tosylchloramide sodium	13103206	8
Caelo	Tosylchloramide sodium	13103207	1
Caelo	Tosylchloramide sodium	71051187	2
Caelo	Tosylchloramide sodium	12323609	2
Caelo	Tosylchloramide sodium	12323607	1
Caelo	Tosylchloramide sodium	12323606	1
Caelo	Tosylchloramide sodium	12323605	1
Caelo	Tosylchloramide sodium	11116708	1
Caelo	Tosylchloramide sodium	14183509	1
Caelo	Tosylchloramide sodium	72215377	2
Caelo	Tragacanth, powdered	9081202	1
Caelo	Tragacanth, powdered	8081201	1
Caelo	Tragacanth, powdered	16110802	2
Caelo	Tragacanth, powdered	16110807	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tragacanth, powdered	1327802	1
Caelo	Tragacanth, powdered	16350004	1
Caelo	Tragacanth, powdered	13277802	4
Caelo	Tragacanth, powdered	14162205	7
Caelo	Tragacanth, powdered	7021402	1
Caelo	Tragacanth, powdered	15085602	4
Caelo	Tragacanth, powdered	153799	1
Caelo	Tragacanth, powdered	20111305	1
Caelo	Tragacanth, powdered	13204701	4
Caelo	Tragacanth, powdered	15085602	1
Caelo	Tragacanth, powdered	12108806	1
Caelo	Tragacanth, powdered	12285408	1
Caelo	Tragacanth, powdered	13095203	1
Caelo	Tragacanth, powdered	15041302	1
Caelo	Tragacanth, powdered	1228503	1
Caelo	Tragacanth, powdered	12285403	3
Caelo	Tragacanth, powdered	14011302	1
Caelo	Tragacanth, powdered	12285402	1
Caelo	Tragacanth, powdered	11128905	1
Caelo	Tragacanth, powdered	15085606	6
Caelo	Tragacanth, powdered	14162207	1
Caelo	Tragacanth, powdered	15121203	1
Caelo	Tragacanth, powdered	15085607	1
Caelo	Tragacanth, powdered	13095208	3
Caelo	Tragacanth, powdered	13277804	3
Caelo	Tretinoin	14040507	1
Caelo	Tretinoin	15239502	1
Caelo	Tretinoin	16217107	1
Caelo	Tretinoin	15078307	3
Caelo	Tretinoin	15H31-B01	1
Caelo	Tretinoin	151000701	1
Caelo	Tretinoin	16291502	5
Caelo	Tretinoin	15437301	1
Caelo	Tretinoin	16F14-B04-324636	2
Caelo	Tretinoin	16070001	1
Caelo	Tretinoin	16291503	6
Caelo	Tretinoin	15078308	1
Caelo	Tretinoin	14252710	2
Caelo	Tretinoin	15239501	2
Caelo	Tretinoin	15239503	2
Caelo	Tretinoin	15078303	1
Caelo	Tretinoin	15078302	1
Caelo	Tretinoin	14252801	1
Caelo	Tretinoin	14252702	3
Caelo	Tretinoin	1612010-01	26
Caelo	Tretinoin	1452802	1
Caelo	Tretinoin	15078302	1
Caelo	Tretinoin	14252706	1
Caelo	Tretinoin	1517E-02563	1
Caelo	Tretinoin	14040504	1
Caelo	Tretinoin	16070003	1
Caelo	Tretinoin	14252706	1
Caelo	Tretinoin	15239502	4
Caelo	Tretinoin	15078301	1
Caelo	Tretinoin	5012I-02563	1
Caelo	Tretinoin	14040504	5
Caelo	Tretinoin	11284703	1
Caelo	Tretinoin	11225506	1
Caelo	Tretinoin	12102904	2
Caelo	Tretinoin	12102902	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tretinoin	12023302	12
Caelo	Tretinoin	171061025343	1
Caelo	Tretinoin	16217120	1
Caelo	Tretinoin	172773102	1
Caelo	Tretinoin	17305701	5
Caelo	Tretinoin	17305702	5
Caelo	Tretinoin	17357001	6
Caelo	Tretinoin	17305707	4
Caelo	Tretinoin	17305704	2
Caelo	Tretinoin	1735707	1
Caelo	Tretinoin	173057075597	1
Caelo	Tretinoin	17357005	2
Caelo	Tretinoin	17305708	1
Caelo	Tretinoin	19050001	2
Caelo	Tretinoin	18167410	1
Caelo	Tretinoin	18167309	7
Caelo	Tretinoin	18167408	1
Caelo	Tretinoin	18167403	1
Caelo	Tretinoin	18167301	1
Caelo	Tretinoin	11225504	1
Caelo	Tretinoin	18167307	1
Caelo	Tretinoin	17305710	1
Caelo	Tretinoin	17357009	2
Caelo	Tretinoin	18167401	5
Caelo	Tretinoin	17357008	3
Caelo	Tretinoin	17305711	1
Caelo	Tretinoin	17305709	1
Caelo	Tretinoin	17357002	1
Caelo	Tretinoin	16K03-806328-388	1
Caelo	Tretinoin	711212563	1
Caelo	Tretinoin	11225511	1
Caelo	Tretinoin	12367104	1
Caelo	Tretinoin	13278101	1
Caelo	Tretinoin	13195002	15
Caelo	Tretinoin	13278104	2
Caelo	Tretinoin	13278102	4
Caelo	Tretinoin	13195003	1
Caelo	Tretinoin	16217108	3
Caelo	Tretinoin	13278105	2
Caelo	Tretinoin	17106102	7
Caelo	Tretinoin	131950002	1
Caelo	Tretinoin	13308003	1
Caelo	Tretinoin	13195008	1
Caelo	Tretinoin	13308004	5
Caelo	Tretinoin	14093502	4
Caelo	Tretinoin	17106101	5
Caelo	Tretinoin	14093503	1
Caelo	Tretinoin	17106103	1
Caelo	Tretinoin	12367101	2
Caelo	Tretinoin	12I06-N01	3
Caelo	Tretinoin	12320201	2
Caelo	Tretinoin	12D05-N08	1
Caelo	Tretinoin	12023309	1
Caelo	Tretinoin	3041302	1
Caelo	Tretinoin	12320205	5
Caelo	Tretinoin	12023303	1
Caelo	Tretinoin	12320301	1
Caelo	Tretinoin	12102903	2
Caelo	Tretinoin	12320203	2
Caelo	Tretinoin	12320207	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Tretinoin	12367103	3
Caelo	Tretinoin	12320209	6
Caelo	Tretinoin	12320303	1
Caelo	Tretinoin	12320204	1
Caelo	Tretinoin	15437303	14
Caelo	Triamcinolone acetonide	18100221	1
Caelo	Triamcinolone acetonide	160102101	1
Caelo	Triamcinolone acetonide	17120611	17
Caelo	Triamcinolone acetonide	3041206	1
Caelo	Triamcinolone acetonide	13124808	12
Caelo	Triamcinolone acetonide	13124808	1
Caelo	Triamcinolone acetonide	1206008-01	1
Caelo	Triamcinolone acetonide	13124811	4
Caelo	Triamcinolone acetonide	13379903	2
Caelo	Triamcinolone acetonide	13256604	6
Caelo	Triamcinolone acetonide	1514Q2516	1
Caelo	Triamcinolone acetonide	17057909	2
Caelo	Triamcinolone acetonide	17120606	14
Caelo	Triamcinolone acetonide	16079502	1
Caelo	Triamcinolone acetonide	17120604	6
Caelo	Triamcinolone acetonide	95/AKD34	1
Caelo	Triamcinolone acetonide	17120605	4
Caelo	Triamcinolone acetonide	13124801	2
Caelo	Triamcinolone acetonide	17057907	4
Caelo	Triamcinolone acetonide	170303701	2
Caelo	Triamcinolone acetonide	6051401	1
Caelo	Triamcinolone acetonide	16K07-B06-333438	6
Caelo	Triamcinolone acetonide	13379913	4
Caelo	Triamcinolone acetonide	13256608	2
Caelo	Triamcinolone acetonide	13379904	4
Caelo	Triamcinolone acetonide	13379914	4
Caelo	Triamcinolone acetonide	13256603	1
Caelo	Triamcinolone acetonide	13379906	5
Caelo	Triamcinolone acetonide	13379912	8
Caelo	Triamcinolone acetonide	13K15-B31-289027	6
Caelo	Triamcinolone acetonide	17120603	13
Caelo	Triamcinolone acetonide	14188302	5
Caelo	Triamcinolone acetonide	13124805	4
Caelo	Triamcinolone acetonide	17120619	6
Caelo	Triamcinolone acetonide	12301508	4
Caelo	Triamcinolone acetonide	120904002	1
Caelo	Triamcinolone acetonide	12301503	2
Caelo	Triamcinolone acetonide	12301504	1
Caelo	Triamcinolone acetonide	12301509	5
Caelo	Triamcinolone acetonide	93/AD1	1
Caelo	Triamcinolone acetonide	1201506	1
Caelo	Triamcinolone acetonide	12301515	9
Caelo	Triamcinolone acetonide	2513M02516	1
Caelo	Triamcinolone acetonide	12301507	4
Caelo	Triamcinolone acetonide	12301518	5
Caelo	Triamcinolone acetonide	13E15-N03	19
Caelo	Triamcinolone acetonide	9081305	1
Caelo	Triamcinolone acetonide	16094507	3
Caelo	Triamcinolone acetonide	13124804	8
Caelo	Triamcinolone acetonide	13e15-n03	1
Caelo	Triamcinolone acetonide	1210-n09	1
Caelo	Triamcinolone acetonide	13124803	2
Caelo	Triamcinolone acetonide	1811211	1
Caelo	Triamcinolone acetonide	4613E2516	1
Caelo	Triamcinolone acetonide	13256607	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triamcinolone acetonide	13124810	1
Caelo	Triamcinolone acetonide	1312407	1
Caelo	Triamcinolone acetonide	13256607	3
Caelo	Triamcinolone acetonide	13124804	3
Caelo	Triamcinolone acetonide	12301515	2
Caelo	Triamcinolone acetonide	13124807	9
Caelo	Triamcinolone acetonide	13124809	1
Caelo	Triamcinolone acetonide	13124806	2
Caelo	Triamcinolone acetonide	13379901	1
Caelo	Triamcinolone acetonide	13379911	1
Caelo	Triamcinolone acetonide	14188304	3
Caelo	Triamcinolone acetonide	15233904	13
Caelo	Triamcinolone acetonide	14355701	1
Caelo	Triamcinolone acetonide	170324T1612014-01	1
Caelo	Triamcinolone acetonide	14233904	1
Caelo	Triamcinolone acetonide	14G23_B07_2976618	1
Caelo	Triamcinolone acetonide	15233905	1
Caelo	Triamcinolone acetonide	16094512	2
Caelo	Triamcinolone acetonide	15233906	1
Caelo	Triamcinolone acetonide	15233903	22
Caelo	Triamcinolone acetonide	310320161	1
Caelo	Triamcinolone acetonide	15233908	2
Caelo	Triamcinolone acetonide	14355712	6
Caelo	Triamcinolone acetonide	16094514	14
Caelo	Triamcinolone acetonide	14355706	1
Caelo	Triamcinolone acetonide	16094504	4
Caelo	Triamcinolone acetonide	15L21-B03-316123	8
Caelo	Triamcinolone acetonide	15430103	7
Caelo	Triamcinolone acetonide	15430105	6
Caelo	Triamcinolone acetonide	5430105	1
Caelo	Triamcinolone acetonide	16094505	2
Caelo	Triamcinolone acetonide	12301507	2
Caelo	Triamcinolone acetonide	160501001	4
Caelo	Triamcinolone acetonide	15430106	5
Caelo	Triamcinolone acetonide	15L21-B03-318663	3
Caelo	Triamcinolone acetonide	16092504	18
Caelo	Triamcinolone acetonide	16094511	15
Caelo	Triamcinolone acetonide	16094506	6
Caelo	Triamcinolone acetonide	160800501	2
Caelo	Triamcinolone acetonide	15233906	16
Caelo	Triamcinolone acetonide	14355102	1
Caelo	Triamcinolone acetonide	1504007-2	1
Caelo	Triamcinolone acetonide	1612014	1
Caelo	Triamcinolone acetonide	13379919	1
Caelo	Triamcinolone acetonide	14355103	13
Caelo	Triamcinolone acetonide	13379921	1
Caelo	Triamcinolone acetonide	14188303	3
Caelo	Triamcinolone acetonide	14355101	12
Caelo	Triamcinolone acetonide	13379922	1
Caelo	Triamcinolone acetonide	14355104	3
Caelo	Triamcinolone acetonide	12301505	1
Caelo	Triamcinolone acetonide	12301502	1
Caelo	Triamcinolone acetonide	14355708	10
Caelo	Triamcinolone acetonide	1355101	1
Caelo	Triamcinolone acetonide	14355702	11
Caelo	Triamcinolone acetonide	16094513	1
Caelo	Triamcinolone acetonide	170547902	1
Caelo	Triamcinolone acetonide	14355711	8
Caelo	Triamcinolone acetonide	14355703	3
Caelo	Triamcinolone acetonide	14355704	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triamcinolone acetonide	17057905	11
Caelo	Triamcinolone acetonide	15233902	21
Caelo	Triamcinolone acetonide	15233901	9
Caelo	Triamcinolone acetonide	17057908	16
Caelo	Triamcinolone acetonide	15C19-B06-308573	10
Caelo	Triamcinolone acetonide	16K07-B06-329447	1
Caelo	Triamcinolone acetonide	17057903	2
Caelo	Triamcinolone acetonide	1508008-01	14
Caelo	Triamcinolone acetonide	1530102	1
Caelo	Triamcinolone acetonide	17057901	7
Caelo	Triamcinolone acetonide	15233905	5
Caelo	Triamcinolone acetonide	17057902	9
Caelo	Triamcinolone acetonide	12301506	2
Caelo	Triamcinolone acetonide	12173704	2
Caelo	Triamcinolone acetonide	16094501	9
Caelo	Triamcinolone acetonide	12173710	4
Caelo	Triamcinolone acetonide	18100226	2
Caelo	Triamcinolone acetonide	9364429	2
Caelo	Triamcinolone acetonide	17120620	4
Caelo	Triamcinolone acetonide	18F14B01355979	1
Caelo	Triamcinolone acetonide	18100220	5
Caelo	Triamcinolone acetonide	18F14-B01-351889	1
Caelo	Triamcinolone acetonide	18100207	1
Caelo	Triamcinolone acetonide	18100213	1
Caelo	Triamcinolone acetonide	18100209	3
Caelo	Triamcinolone acetonide	18100225	1
Caelo	Triamcinolone acetonide	42(18100213)	1
Caelo	Triamcinolone acetonide	18100203	3
Caelo	Triamcinolone acetonide	18100205	2
Caelo	Triamcinolone acetonide	17358302	14
Caelo	Triamcinolone acetonide	18100204	5
Caelo	Triamcinolone acetonide	17358302-T20181115Ma	1
Caelo	Triamcinolone acetonide	17358303	4
Caelo	Triamcinolone acetonide	18100206	4
Caelo	Triamcinolone acetonide	18010202-01	1
Caelo	Triamcinolone acetonide	1720608	1
Caelo	Triamcinolone acetonide	18100210	8
Caelo	Triamcinolone acetonide	17356306	1
Caelo	Triamcinolone acetonide	17358304	7
Caelo	Triamcinolone acetonide	17128606	1
Caelo	Triamcinolone acetonide	17120610	8
Caelo	Triamcinolone acetonide	12173707	5
Caelo	Triamcinolone acetonide	17120608	8
Caelo	Triamcinolone acetonide	11267720	1
Caelo	Triamcinolone acetonide	11267704	1
Caelo	Triamcinolone acetonide	10220341	1
Caelo	Triamcinolone acetonide	12173705	2
Caelo	Triamcinolone acetonide	1107212516	1
Caelo	Triamcinolone acetonide	11267710	2
Caelo	Triamcinolone acetonide	17120618	6
Caelo	Triamcinolone acetonide	1706601-01	1
Caelo	Triamcinolone acetonide	NS0571M	6
Caelo	Triamcinolone acetonide	11267721	2
Caelo	Triamcinolone acetonide	11267707	4
Caelo	Triamcinolone acetonide	12173701	5
Caelo	Triamcinolone, micron...	15106902	1
Caelo	Triamcinolone, micron...	12295210	1
Caelo	Triamcinolone, micron...	14B13-B30-293474	2
Caelo	Triamcinolone, micron...	15106960	1
Caelo	Triamcinolone, micron...	15106912	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triamcinolone, micron...	13J07-N08	2
Caelo	Triamcinolone, micron...	18092104	1
Caelo	Trichloroacetic acid	12135203	1
Caelo	Trichloroacetic acid	13056213	1
Caelo	Trichloroacetic acid	13056212	1
Caelo	Trichloroacetic acid	6091306	1
Caelo	Trichloroacetic acid	13056209	1
Caelo	Trichloroacetic acid	12135210	2
Caelo	Trichloroacetic acid	10379401	1
Caelo	Trichloroacetic acid	13056203	3
Caelo	Trichloroacetic acid	12135207	5
Caelo	Trichloroacetic acid	11294318	1
Caelo	Trichloroacetic acid	12g04m11	1
Caelo	Trichloroacetic acid	14101711	1
Caelo	Trichloroacetic acid	12135201	1
Caelo	Trichloroacetic acid	10130802	1
Caelo	Triclosan	1609005	1
Caelo	Triclosan	18081408	15
Caelo	Triclosan	16084618	14
Caelo	Triclosan	17149813-T20181115Ma	1
Caelo	Triclosan	18081407	11
Caelo	Triclosan	16084614	19
Caelo	Triclosan	18081406	9
Caelo	Triclosan	16084617	1
Caelo	Triclosan	18253706	8
Caelo	Triclosan	10206240	1
Caelo	Triclosan	1714926	1
Caelo	Triclosan	16084608	17
Caelo	Triclosan	16084606	4
Caelo	Triclosan	18253707	8
Caelo	Triclosan	17149822	8
Caelo	Triclosan	1719814	1
Caelo	Triclosan	2317I-03402	1
Caelo	Triclosan	18081410	5
Caelo	Triclosan	17081709	13
Caelo	Triclosan	17149819	11
Caelo	Triclosan	18081414	9
Caelo	Triclosan	1801417	1
Caelo	Triclosan	16084611	15
Caelo	Triclosan	18081417	13
Caelo	Triclosan	18253701	2
Caelo	Triclosan	162820817	1
Caelo	Triclosan	1808414	1
Caelo	Triclosan	15C10-B06-310031	1
Caelo	Triclosan	18081413	16
Caelo	Triclosan	16280817	12
Caelo	Triclosan	18081415	11
Caelo	Triclosan	17081711	14
Caelo	Triclosan	12276827	1
Caelo	Triclosan	12145909	10
Caelo	Triclosan	17149815	9
Caelo	Triclosan	11238003	2
Caelo	Triclosan	18081403	8
Caelo	Triclosan	1714930	1
Caelo	Triclosan	15301513	1
Caelo	Triclosan	160084601	1
Caelo	Triclosan	16084601	6
Caelo	Triclosan	15301511	7
Caelo	Triclosan	1530105	1
Caelo	Triclosan	16084607	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triclosan	1503603	1
Caelo	Triclosan	15301504	10
Caelo	Triclosan	17149820	13
Caelo	Triclosan	171498085598	1
Caelo	Triclosan	15301512	5
Caelo	Triclosan	1714819	1
Caelo	Triclosan	17149824	5
Caelo	Triclosan	17149825	9
Caelo	Triclosan	17K02-B06_343544	2
Caelo	Triclosan	15301514	1
Caelo	Triclosan	16084602	8
Caelo	Triclosan	18081401	6
Caelo	Triclosan	18253702	8
Caelo	Triclosan	16084609	3
Caelo	Triclosan	17149814	9
Caelo	Triclosan	16071801	1
Caelo	Triclosan	20021304	1
Caelo	Triclosan	18081402	6
Caelo	Triclosan	17149827	8
Caelo	Triclosan	160846	1
Caelo	Triclosan	16084610	8
Caelo	Triclosan	1748912	1
Caelo	Triclosan	16084603	10
Caelo	Triclosan	17149812-20180619Ma4	1
Caelo	Triclosan	170802-01	1
Caelo	Triclosan	17149808	21
Caelo	Triclosan	17149821	17
Caelo	Triclosan	182537052	1
Caelo	Triclosan	160846018	1
Caelo	Triclosan	18253705	3
Caelo	Triclosan	17081722	22
Caelo	Triclosan	228711-514	1
Caelo	Triclosan	6121711	1
Caelo	Triclosan	17149806	8
Caelo	Triclosan	1149806	1
Caelo	Triclosan	12003205	7
Caelo	Triclosan	17081724	15
Caelo	Triclosan	1628082	1
Caelo	Triclosan	1701032-01	20
Caelo	Triclosan	16280803	4
Caelo	Triclosan	16280816	12
Caelo	Triclosan	16280819	7
Caelo	Triclosan	16280820	9
Caelo	Triclosan	17081713	4
Caelo	Triclosan	16280824	2
Caelo	Triclosan	17081705	8
Caelo	Triclosan	16280874	1
Caelo	Triclosan	17081725	14
Caelo	Triclosan	11076716	1
Caelo	Triclosan	17149805	5
Caelo	Triclosan	11238009	2
Caelo	Triclosan	12c21-n05	4
Caelo	Triclosan	12003201	5
Caelo	Triclosan	170817075274	1
Caelo	Triclosan	17081712	12
Caelo	Triclosan	17081710	15
Caelo	Triclosan	17081707	24
Caelo	Triclosan	170501101	3
Caelo	Triclosan	16280818	7
Caelo	Triclosan	17081706	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triclosan	17081717	3
Caelo	Triclosan	16280821	2
Caelo	Triclosan	17081704	1
Caelo	Triclosan	16280809	5
Caelo	Triclosan	18253709	5
Caelo	Triclosan	11076717	2
Caelo	Triclosan	17149801	6
Caelo	Triclosan	12145910	13
Caelo	Triclosan	12145908	8
Caelo	Triclosan	170802501	2
Caelo	Triclosan	170817265507	1
Caelo	Triclosan	17149812	8
Caelo	Triclosan	18253716	1
Caelo	Triclosan	16280810	1
Caelo	Triclosan	17149817	8
Caelo	Triclosan	12680801	1
Caelo	Triclosan	17081721	13
Caelo	Triclosan	18253712	3
Caelo	Triclosan	18253715	1
Caelo	Triclosan	16280801	25
Caelo	Triclosan	162808	1
Caelo	Triclosan	16084620	1
Caelo	Triclosan	17149809	25
Caelo	Triclosan	17081720	14
Caelo	Triclosan	16280808	17
Caelo	Triclosan	17081726	7
Caelo	Triclosan	17081716	12
Caelo	Triclosan	16280805	1
Caelo	Triclosan	17081716-T20180112Ma	1
Caelo	Triclosan	17081714	7
Caelo	Triclosan	16280811	16
Caelo	Triclosan	12145901	5
Caelo	Triclosan	12145902	5
Caelo	Triclosan	5012E-03402	1
Caelo	Triclosan	17149813	13
Caelo	Triclosan	16280802	18
Caelo	Triclosan	17149807	11
Caelo	Triclosan	170817	1
Caelo	Triclosan	12145906	6
Caelo	Triclosan	17081703	10
Caelo	Triclosan	12003210	5
Caelo	Triclosan	517E-03402	1
Caelo	Triclosan	1708173	1
Caelo	Triclosan	ch.b.15301506	1
Caelo	Triclosan	13D02-N02	80
Caelo	Triclosan	14211207	1
Caelo	Triclosan	1412109	1
Caelo	Triclosan	13L18-B35-289992	5
Caelo	Triclosan	14121209	8
Caelo	Triclosan	14121203	8
Caelo	Triclosan	15203608	11
Caelo	Triclosan	1520603	1
Caelo	Triclosan	15203604	12
Caelo	Triclosan	14052105	2
Caelo	Triclosan	15G07-B03-30039	1
Caelo	Triclosan	15203607	16
Caelo	Triclosan	15301501	1
Caelo	Triclosan	14121207	11
Caelo	Triclosan	15301502	14
Caelo	Triclosan	15G07-B03-310035	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Triclosan	5203607	1
Caelo	Triclosan	12359110	5
Caelo	Triclosan	12359104	8
Caelo	Triclosan	14121204	9
Caelo	Triclosan	13L18-B35-289995	12
Caelo	Triclosan	23121409	1
Caelo	Triclosan	14121215	16
Caelo	Triclosan	12359108	13
Caelo	Triclosan	14205901	5
Caelo	Triclosan	12359111	9
Caelo	Triclosan	14205904	7
Caelo	Triclosan	14205902	9
Caelo	Triclosan	150103201	1
Caelo	Triclosan	12359105	10
Caelo	Triclosan	14205905	5
Caelo	Triclosan	1505009	1
Caelo	Triclosan	15203603	17
Caelo	Triclosan	150704101	1
Caelo	Triclosan	14205903	12
Caelo	Triclosan	12359106	13
Caelo	Triclosan	14211211	1
Caelo	Triclosan	12359112	9
Caelo	Triclosan	11777748	7
Caelo	Triclosan	14121214	10
Caelo	Triclosan	14121211	10
Caelo	Triclosan	13L18B35289986	1
Caelo	Triclosan	14121220	1
Caelo	Triclosan	8081401	1
Caelo	Triclosan	14F23B02293859	1
Caelo	Triclosan	15301505	11
Caelo	Triclosan	13247805	1
Caelo	Triclosan	13D02-No2	1
Caelo	Triclosan	13391801	2
Caelo	Triclosan	15121202	1
Caelo	Triclosan	12359109	4
Caelo	Triclosan	AKD12	1
Caelo	Triclosan	12145907	9
Caelo	Triclosan	3041408	1
Caelo	Triclosan	13247905	11
Caelo	Triclosan	16021306	1
Caelo	Triclosan	13391802	4
Caelo	Triclosan	1125204	1
Caelo	Triclosan	11238011	1
Caelo	Triclosan	14052103	14
Caelo	Triclosan	13247908	2
Caelo	Triclosan	13L18B35289993	1
Caelo	Triclosan	12145917	1
Caelo	Triclosan	15301503	9
Caelo	Triclosan	12L06-N01/0005900346	1
Caelo	Triclosan	15301506	15
Caelo	Triclosan	15G07-B03-310035/201	1
Caelo	Triclosan	14052102	12
Caelo	Triclosan	13391808	9
Caelo	Triclosan	15301510	3
Caelo	Triclosan	14052101	5
Caelo	Triclosan	bhfdd	1
Caelo	Trometamol	16168820	16
Caelo	Trometamol	17197804	14
Caelo	Trometamol	16168801	1
Caelo	Trometamol	1618810	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Trometamol	16D14-b-09-328729	1
Caelo	Trometamol	13095913	3
Caelo	Trometamol	17053184	1
Caelo	Trometamol	314Q2519	1
Caelo	Trometamol	16168825	2
Caelo	Trometamol	13095912	17
Caelo	Trometamol	17320801	20
Caelo	Trometamol	27111304	1
Caelo	Trometamol	17320816	8
Caelo	Trometamol	15377602	19
Caelo	Trometamol	18197101	4
Caelo	Trometamol	15277602	1
Caelo	Trometamol	15377607	28
Caelo	Trometamol	1730804	1
Caelo	Trometamol	17197806	2
Caelo	Trometamol	17320801/1160518	1
Caelo	Trometamol	18041815	1
Caelo	Trometamol	17197809	5
Caelo	Trometamol	171978	1
Caelo	Trometamol	181016182	1
Caelo	Trometamol	16168810	29
Caelo	Trometamol	17020304	7
Caelo	Trometamol	15070309	1
Caelo	Trometamol	150703	1
Caelo	Trometamol	16168812	2
Caelo	Trometamol	15070307	4
Caelo	Trometamol	15377611	22
Caelo	Trometamol	215I-02519	1
Caelo	Trometamol	15070308	24
Caelo	Trometamol	15094404	1
Caelo	Trometamol	14050003	3
Caelo	Trometamol	15070304	9
Caelo	Trometamol	18197103	4
Caelo	Trometamol	14050004	3
Caelo	Trometamol	18081503	1
Caelo	Trometamol	14050006	1
Caelo	Trometamol	15377605	1
Caelo	Trometamol	14332404	1
Caelo	Trometamol	914I-02519	1
Caelo	Trometamol	14332402	23
Caelo	Trometamol	15094401	13
Caelo	Trometamol	23062015B	1
Caelo	Trometamol	3713E-02519	1
Caelo	Trometamol	16168807	1
Caelo	Trometamol	13C-08-N07	1
Caelo	Trometamol	13339108	11
Caelo	Trometamol	5315I-02519	1
Caelo	Trometamol	ch.b.15377602	1
Caelo	Trometamol	4616M-02519	1
Caelo	Trometamol	161013T/15377607	1
Caelo	Trometamol	8031301	1
Caelo	Trometamol	13095915	4
Caelo	Trometamol	1180116	1
Caelo	Trometamol	15070311	21
Caelo	Trometamol	14332403	1
Caelo	Trometamol	13339103	7
Caelo	Trometamol	5012M-02519	1
Caelo	Trometamol	13339107	7
Caelo	Trometamol	1339103	1
Caelo	Trometamol	13051401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Trometamol	15377604	2
Caelo	Trometamol	14050002	7
Caelo	Trometamol	17320810	1
Caelo	Trometamol	16168816	2
Caelo	Trometamol	12113510	2
Caelo	Trometamol	12113511	1
Caelo	Trometamol	11376006	8
Caelo	Trometamol	26111201	1
Caelo	Trometamol	11160707	1
Caelo	Trometamol	213E-02519	1
Caelo	Trometamol	1001312519	1
Caelo	Trometamol	13095911	13
Caelo	Trometamol	11376005	6
Caelo	Trometamol	13095905	1
Caelo	Trometamol	13095910	6
Caelo	Trometamol	12113507	3
Caelo	Trometamol	11160711	1
Caelo	Trometamol	11160706	1
Caelo	Trometamol	12113502	3
Caelo	Trometamol	12113505	1
Caelo	Trometamol		2
Caelo	Trometamol	21082013A	1
Caelo	Trometamol	12113515	4
Caelo	Trometamol	17320611	1
Caelo	Trometamol	13095902	13
Caelo	Trometamol	13095901	9
Caelo	Tryptophan (L-) (H1053)	15000441/4	2
Caelo	Turmeric root, cut	12136506	1
Caelo	Turmeric root, cut	12136507	1
Caelo	Turmeric root, cut	15158001	2
Caelo	Turmeric root, cut	13360204	2
Caelo	Turmeric root, cut	15372603	1
Caelo	Turmeric root, cut	1213650	1
Caelo	Turmeric root, cut	15207601	1
Caelo	Turmeric root, cut	12136503	2
Caelo	Turmeric root, cut	12136504	1
Caelo	Turmeric root, cut	12136501	2
Caelo	Turmeric root, cut	15158001	2
Caelo	Turmeric, javanese, cut	17042003	1
Caelo	Turmeric, javanese, cut	13262907	1
Caelo	Turmeric, javanese, cut	15120202	2
Caelo	Turmeric, javanese, cut	17042002	1
Caelo	Turmeric, javanese, cut	12319206	2
Caelo	Turmeric, javanese, cut	12319201	1
Caelo	Turmeric, javanese, cut	12319206	2
Caelo	Turmeric, javanese, cut	14309503	1
Caelo	Turmeric, javanese, cut	25041313	1
Caelo	Turmeric, javanese, cut	11182903	1
Caelo	Turmeric, javanese, cut	14017706	1
Caelo	Ubidecarenone	18048401	1
Caelo	Urea	16107425	5
Caelo	Urea	14107428	1
Caelo	Urea	16107427	1
Caelo	Urea	18024201	26
Caelo	Urea	1802402	1
Caelo	Urea	16107429	6
Caelo	Urea	16107423	11
Caelo	Urea	21041802	1
Caelo	Urea	1607422	1
Caelo	Urea	180322H/18024202	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Urea	12L10-N08	6
Caelo	Urea	1241017	1
Caelo	Urea	12111803	1
Caelo	Urea	18024202, 180425UP	1
Caelo	Urea	16107428	22
Caelo	Urea	2324534	1
Caelo	Urea	12111834	2
Caelo	Urea	26071304	1
Caelo	Urea	13064201	7
Caelo	Urea	107063	1
Caelo	Urea	13064206	22
Caelo	Urea	16107248	1
Caelo	Urea	13064203	13
Caelo	Urea	123456	2
Caelo	Urea	1211820	1
Caelo	Urea	21051304	1
Caelo	Urea	12111826	7
Caelo	Urea	1211809	1
Caelo	Urea	16107426	9
Caelo	Urea	17107429	1
Caelo	Urea	3012A-02544	1
Caelo	Urea	11150714	2
Caelo	Urea	21111806	1
Caelo	Urea	12111820	13
Caelo	Urea	18023010	17
Caelo	Urea	20130912-4	1
Caelo	Urea	18109801	1
Caelo	Urea	18023011	28
Caelo	Urea	18023012	21
Caelo	Urea	12111809	12
Caelo	Urea	84(18023010)	1
Caelo	Urea	1211808	1
Caelo	Urea	12111807	5
Caelo	Urea	18023015	14
Caelo	Urea	85(18023012)	2
Caelo	Urea	10169726	2
Caelo	Urea	18203012	1
Caelo	Urea	11150702	2
Caelo	Urea	18023016	2
Caelo	Urea	4111M-02544	1
Caelo	Urea	18023018	1
Caelo	Urea	10116	1
Caelo	Urea	19031604	2
Caelo	Urea	1211804	1
Caelo	Urea	12111822	20
Caelo	Urea	1211806	1
Caelo	Urea	17321206	1
Caelo	Urea	18023003-20180830Ma2	1
Caelo	Urea	18023001	9
Caelo	Urea	30031304	1
Caelo	Urea	1802401	1
Caelo	Urea	16/13-B02-332593	1
Caelo	Urea	18023003	23
Caelo	Urea	13064226	12
Caelo	Urea	1061801	1
Caelo	Urea	6101703	1
Caelo	Urea	10169722	1
Caelo	Urea	18023007	16
Caelo	Urea	18024202	29
Caelo	Urea	18023004	66

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Urea	1803003	1
Caelo	Urea	11150725	3
Caelo	Urea	1823004	1
Caelo	Urea	12111802	8
Caelo	Urea	180223003	1
Caelo	Urea	10081208	1
Caelo	Urea	118023004	1
Caelo	Urea	2109212544	1
Caelo	Urea	18023006	31
Caelo	Urea	12111806	8
Caelo	Urea	12111823	12
Caelo	Urea	12111801	5
Caelo	Urea	11150722	3
Caelo	Urea	1364212	1
Caelo	Urea	Val. Mes.: 4415I02544	1
Caelo	Urea	13064207	18
Caelo	Urea	1514110	2
Caelo	Urea	14107120	8
Caelo	Urea	14107118	7
Caelo	Urea	160620UP	1
Caelo	Urea	4116M-02544	1
Caelo	Urea	15154116	13
Caelo	Urea	1236	1
Caelo	Urea	1514107	1
Caelo	Urea	16107407	29
Caelo	Urea	13042015B	1
Caelo	Urea	151541111	1
Caelo	Urea	16107405	8
Caelo	Urea	14107119	21
Caelo	Urea	1610407	1
Caelo	Urea	10038587	2
Caelo	Urea	16107412	75
Caelo	Urea	1514112	1
Caelo	Urea	16107413	8
Caelo	Urea	161074	2
Caelo	Urea	16107411	33
Caelo	Urea	1517E-02544	1
Caelo	Urea	63/AKD4	1
Caelo	Urea	14107122	15
Caelo	Urea	15087102	7
Caelo	Urea	15087103	22
Caelo	Urea	15154103	49
Caelo	Urea	17.02.16Vi	1
Caelo	Urea	Validierungsmessung	1
Caelo	Urea	15154107	24
Caelo	Urea	160125/2	1
Caelo	Urea	15154101	1
Caelo	Urea	151541	2
Caelo	Urea	2/20141507	1
Caelo	Urea	15154109	7
Caelo	Urea	15154111	43
Caelo	Urea	19041703	1
Caelo	Urea	4914I-02544	1
Caelo	Urea	2416I-02544	1
Caelo	Urea	15154112	38
Caelo	Urea	1508/104	1
Caelo	Urea	15087104	17
Caelo	Urea	14197118	1
Caelo	Urea	568900	1
Caelo	Urea	14/18-B03-303581	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Urea	15087107	18
Caelo	Urea	25082015A	1
Caelo	Urea	13072015B	1
Caelo	Urea	2116A-02544	1
Caelo	Urea	20031514	1
Caelo	Urea	14107117	11
Caelo	Urea	14107108	12
Caelo	Urea	13064232	15
Caelo	Urea	1113A-025	1
Caelo	Urea	16107419	19
Caelo	Urea	16107420	36
Caelo	Urea	16107422	27
Caelo	Urea	13064228	22
Caelo	Urea	62/AKD49	1
Caelo	Urea	13064221	6
Caelo	Urea	314Q2544	1
Caelo	Urea	17011405	1
Caelo	Urea	13064227	6
Caelo	Urea	13064223	36
Caelo	Urea	13064224	2
Caelo	Urea	13064213	5
Caelo	Urea	13064215	23
Caelo	Urea	9240908	1
Caelo	Urea	13E21-N05	3
Caelo	Urea	1306412	1
Caelo	Urea	17101310	1
Caelo	Urea	AKD6	1
Caelo	Urea	11150703	1
Caelo	Urea	161072422	1
Caelo	Urea	13064214	13
Caelo	Urea	4013Q2544	1
Caelo	Urea	1306422	1
Caelo	Urea	10038593	1
Caelo	Urea	14107112	20
Caelo	Urea	14011513	1
Caelo	Urea	5208	1
Caelo	Urea	14107113	19
Caelo	Urea	14107109	22
Caelo	Urea	17021601	4
Caelo	Urea	5100	1
Caelo	Urea	1222	1
Caelo	Urea	16107417	25
Caelo	Urea	16107411,20170427Ma8	1
Caelo	Urea	161074165325	1
Caelo	Urea	14107103	8
Caelo	Urea	14E07-B05-289706	2
Caelo	Urea	14107102	28
Caelo	Urea	13064233	16
Caelo	Urea	170509UP/16107407	1
Caelo	Urea	13064236	9
Caelo	Urea	13064231	12
Caelo	Urea	13064230	15
Caelo	Urea	16107416	19
Caelo	Urea	16107418	27
Caelo	Urea	5061410	1
Caelo	Urea	1317B02-287602	1
Caelo	Urea	12111810	17
Caelo	Urea	16107406	19
Caelo	Valerian root, cut	15175205	1
Caelo	Valerian root, cut	11103608	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Valerian root, cut	13287402	3
Caelo	Valerian root, cut	28101303	1
Caelo	Valerian root, cut	11103602	1
Caelo	Valerian root, cut	11103601	1
Caelo	Valerian root, cut	13196801-3564412081	1
Caelo	Valerian root, cut	11103609	4
Caelo	Valerian root, cut	20081303	1
Caelo	Valerian root, cut	13119001	1
Caelo	Valerian root, cut	13196806	6
Caelo	Valerian root, cut	120338012	1
Caelo	Valerian root, cut	13196802	4
Caelo	Valerian root, cut	29101303	1
Caelo	Valerian root, cut	1319003	1
Caelo	Valerian root, cut	12033810	1
Caelo	Valerian root, cut	13119003	12
Caelo	Valerian root, cut	1203381	1
Caelo	Valerian root, cut	11103612	2
Caelo	Valerian root, cut	13119002	18
Caelo	Valerian root, cut	12033811	5
Caelo	Valerian root, cut	12033812	8
Caelo	Valerian root, cut	12033807	1
Caelo	Valerian root, cut	12033809	9
Caelo	Valerian root, cut	11349101	4
Caelo	Valerian root, cut	12033808	7
Caelo	Valerian root, cut	13119004	3
Caelo	Valerian root, cut	2712A-02167	1
Caelo	Valerian root, cut	12033805	4
Caelo	Valerian root, cut	31121201	1
Caelo	Valerian root, cut	11103607	3
Caelo	Valerian root, cut	31101202	1
Caelo	Valerian root, cut	11035002	1
Caelo	Valerian root, cut	3091202	1
Caelo	Valerian root, powdered	13287402	2
Caelo	Valerian root, powdered	15175204	1
Caelo	Vanillin	12079407	2
Caelo	Vanillin	10066811	1
Caelo	Vanillin	12079410	1
Caelo	Vanillin	12079424	2
Caelo	Vanillin	12079417	1
Caelo	Vanillin	90961269	1
Caelo	Vanillin	112901	1
Caelo	Vanillin	15154020	6
Caelo	Vanillin	12079428	1
Caelo	Vanillin	18073507	2
Caelo	Vanillin	18073506	2
Caelo	Vanillin	1089501	1
Caelo	Vanillin	16089508	5
Caelo	Vanillin	12071401	1
Caelo	Vanillin	15227317	1
Caelo	Vanillin	16089506	8
Caelo	Vanillin	14044419	2
Caelo	Vanillin	15154017	1
Caelo	Vanillin	15154006	2
Caelo	Vanillin	15154014	3
Caelo	Vanillin	15135302, PZN10999111	1
Caelo	Vanillin	16089501	3
Caelo	Vanillin	15135302	1
Caelo	Verbena herb, cut	11031105	2
Caelo	Verbena herb, cut	11031106	5
Caelo	Verbena herb, cut	12188202	14

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Verbena herb, cut	11031108	5
Caelo	Verbena herb, cut	11031103	1
Caelo	Verbena herb, cut	11031110	6
Caelo	Verbena herb, cut	19111205	1
Caelo	Verbena herb, cut	110031110	1
Caelo	Verbena herb, cut	12187601	15
Caelo	Verbena herb, cut	12383005	1
Caelo	Verbena herb, cut	12178601	1
Caelo	Verbena herb, cut	12188203	6
Caelo	Verbena herb, cut	4651b120611	1
Caelo	Verbena herb, cut	13151001	17
Caelo	Verbena herb, cut	12383002	2
Caelo	Verbena herb, cut	6111301	1
Caelo	Verbena herb, cut	13151003	5
Caelo	Verbena herb, cut	2x12188202	1
Caelo	Verbena herb, cut	13351401	1
Caelo	Verbena herb, cut	15012405	1
Caelo	Walnut bark, cut	11325606	1
Caelo	Walnut bark, cut	11325603	1
Caelo	Walnut bark, cut	11016605	1
Caelo	Walnut bark, cut	11016604	1
Caelo	Walnut leaf, cut	12124605	6
Caelo	Walnut leaf, cut	14064910	1
Caelo	Walnut leaf, cut	14064909	1
Caelo	Walnut leaf, cut	14064906	3
Caelo	Walnut leaf, cut	14064904	1
Caelo	Walnut leaf, cut	14064903	1
Caelo	Walnut leaf, cut	15126004	1
Caelo	Walnut leaf, cut	12124502	1
Caelo	Walnut leaf, cut	14064902	4
Caelo	Walnut leaf, cut	1414Q01499	1
Caelo	Walnut leaf, cut	414Q1499	1
Caelo	Walnut leaf, cut	11103010	5
Caelo	Walnut leaf, cut	4713U01499	1
Caelo	Walnut leaf, cut	12124603-341930192	1
Caelo	Walnut leaf, cut	1214603	1
Caelo	Walnut leaf, cut	12124601	5
Caelo	Walnut leaf, cut	3051301	1
Caelo	Walnut leaf, cut	11103008	4
Caelo	Walnut leaf, cut	11103007	5
Caelo	Walnut leaf, cut	11103004	2
Caelo	Walnut leaf, cut	11103003	2
Caelo	Walnut leaf, cut	10050416	1
Caelo	Walnut leaf, cut	15126003	6
Caelo	Walnut leaf, cut	12124602	3
Caelo	Walnut leaf, cut	15126008	5
Caelo	Walnut leaf, cut	3712Q-01499	1
Caelo	Wax, bleached	16163114	15
Caelo	Wax, bleached	151210CA	1
Caelo	Wax, bleached	14336616	1
Caelo	Wax, bleached	15292014	19
Caelo	Wax, bleached	15292010	1
Caelo	Wax, bleached	15317609	1
Caelo	Wax, bleached	15292009	3
Caelo	Wax, bleached	15317608	1
Caelo	Wax, bleached	14336621	1
Caelo	Wax, bleached	1529003	1
Caelo	Wax, bleached	15292004	6
Caelo	Wax, bleached	15292003	8
Caelo	Wax, bleached	14336620	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wax, bleached	18048906	14
Caelo	Wax, bleached	18048907	3
Caelo	Wax, bleached	15292005	2
Caelo	Wax, bleached	15292014, 20161112Ma8	1
Caelo	Wax, bleached	16163101	1
Caelo	Wax, bleached	15292002	1
Caelo	Wax, bleached	3051808	1
Caelo	Wax, bleached	143366201.NP(ca50	1
Caelo	Wax, bleached	17091802	18
Caelo	Wax, bleached	18048901	17
Caelo	Wax, bleached	16163115	2
Caelo	Wax, bleached	122570	1
Caelo	Wax, bleached	29091805	1
Caelo	Wax, bleached	16163106	1
Caelo	Wax, bleached	18048903	2
Caelo	Wax, bleached	13378640	2
Caelo	Wax, bleached	16163109	18
Caelo	Wax, bleached	2271200	1
Caelo	Wax, bleached	15292011	1
Caelo	Wax, bleached	16163105	26
Caelo	Wax, bleached	15292015	1
Caelo	Wax, bleached	17091803	5
Caelo	Wax, bleached	18239106	4
Caelo	Wax, bleached	18239102	10
Caelo	Wax, yellow	16038201	19
Caelo	Wax, yellow	11285001	1
Caelo	Wax, yellow	15317602	4
Caelo	Wax, yellow	14336616	9
Caelo	Wax, yellow	15317608	8
Caelo	Wax, yellow	12124810	5
Caelo	Wax, yellow	13378636	1
Caelo	Wax, yellow	12041303	1
Caelo	Wax, yellow	15292014	1
Caelo	Wax, yellow	20031517	1
Caelo	Wax, yellow	15041301	1
Caelo	Wax, yellow	12124814	2
Caelo	Wax, yellow	12257012	4
Caelo	Wax, yellow	12124816	1
Caelo	Wax, yellow	13031201	4
Caelo	Wax, yellow	12367901	14
Caelo	Wax, yellow	1214816	1
Caelo	Wax, yellow	12257004	2
Caelo	Wax, yellow	13378640	1
Caelo	Wax, yellow	1436607	1
Caelo	Wax, yellow	18121201	1
Caelo	Wax, yellow	4910A-01239	1
Caelo	Wax, yellow	13031206	9
Caelo	Wax, yellow	15317601	5
Caelo	Wax, yellow	12124805	4
Caelo	Wax, yellow	14189914	4
Caelo	Wax, yellow	12367911	2
Caelo	Wax, yellow	15022904	4
Caelo	Wax, yellow	15317609	10
Caelo	Wax, yellow	14336609	2
Caelo	Wax, yellow	15317603	4
Caelo	Wax, yellow	15292007	1
Caelo	Wax, yellow	14189907	2
Caelo	Wax, yellow	14189913	2
Caelo	Wax, yellow	26091302	1
Caelo	Wax, yellow	12367910	15

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wax, yellow	10364810	1
Caelo	Wax, yellow	15317610	6
Caelo	Wax, yellow	12257002	3
Caelo	Wax, yellow	12124806	2
Caelo	Wax, yellow	12124812	8
Caelo	Wax, yellow	12367912	3
Caelo	Wax, yellow	15292009	1
Caelo	Wax, yellow	17121201	1
Caelo	Wax, yellow	1436616	1
Caelo	Wax, yellow	14189902	2
Caelo	Wax, yellow	15023003	4
Caelo	Wax, yellow	15317611	1
Caelo	Wax, yellow	150622CF	1
Caelo	Wax, yellow	12367904	6
Caelo	Wax, yellow	14336605	2
Caelo	Wax, yellow	15022908	10
Caelo	Wax, yellow	11285015	1
Caelo	Wax, yellow	12257022	3
Caelo	Wax, yellow	20031516	1
Caelo	Wax, yellow	14189919	1
Caelo	Wax, yellow	11308311	2
Caelo	Wax, yellow	121248024-2015	1
Caelo	Wax, yellow	22071301	1
Caelo	Wax, yellow	12124802	2
Caelo	Wax, yellow	15317613	9
Caelo	Wax, yellow	13061503	1
Caelo	Wax, yellow	14336608	3
Caelo	Wax, yellow	12124801	4
Caelo	Wax, yellow	13335517	4
Caelo	Wax, yellow	14189915	2
Caelo	Wax, yellow	14189912	4
Caelo	Wax, yellow	14189920	1
Caelo	Wax, yellow	17121202	1
Caelo	Wax, yellow	1436620	1
Caelo	Wax, yellow	14189909	5
Caelo	Wax, yellow	15098902	1
Caelo	Wax, yellow	13335516	1
Caelo	Wax, yellow	14189908	1
Caelo	Wax, yellow	13378639	3
Caelo	Wax, yellow	28011501	1
Caelo	Wax, yellow	15292003	2
Caelo	Wax, yellow	161202CF/15317611	1
Caelo	Wax, yellow	1418915	1
Caelo	Wax, yellow	12124817	2
Caelo	Wax, yellow	11285014	2
Caelo	Wax, yellow	10364837	1
Caelo	Wax, yellow	14336620	3
Caelo	Wax, yellow	14336617	1
Caelo	Wax, yellow	14189917	3
Caelo	Wax, yellow	11308305	1
Caelo	Wax, yellow	15317607	5
Caelo	Wax, yellow	12367916-3429311150	1
Caelo	Wax, yellow	30101308	1
Caelo	Wax, yellow	12124807	3
Caelo	Wax, yellow	13378634	6
Caelo	Wax, yellow	5121707	1
Caelo	Wax, yellow	152920004	1
Caelo	Wax, yellow	14336607	11
Caelo	Wax, yellow	13378602	1
Caelo	Wax, yellow	13051401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wax, yellow	15292004	7
Caelo	Wax, yellow	17011611	5
Caelo	Wax, yellow	170011607	2
Caelo	Wax, yellow	11017315	1
Caelo	Wax, yellow	13335507	6
Caelo	Wax, yellow	12367916	6
Caelo	Wax, yellow	1701611	1
Caelo	Wax, yellow	17207402	14
Caelo	Wax, yellow	1595	1
Caelo	Wax, yellow	19041704	1
Caelo	Wax, yellow	12367020	1
Caelo	Wax, yellow	13335508	1
Caelo	Wax, yellow	17041402	1
Caelo	Wax, yellow	13031222	3
Caelo	Wax, yellow	16163112	1
Caelo	Wax, yellow	17011619	1
Caelo	Wax, yellow	17011604	12
Caelo	Wax, yellow	13335509	2
Caelo	Wax, yellow	16038210	6
Caelo	Wax, yellow	12K20B05287421	1
Caelo	Wax, yellow	17011605	1
Caelo	Wax, yellow	5022904	1
Caelo	Wax, yellow	14G07-B04-300495	3
Caelo	Wax, yellow	170011602	1
Caelo	Wax, yellow	16051702	1
Caelo	Wax, yellow	17011602	5
Caelo	Wax, yellow	14336618	1
Caelo	Wax, yellow	17011609	10
Caelo	Wax, yellow	17011607	9
Caelo	Wax, yellow	20071503	1
Caelo	Wax, yellow	1701607	1
Caelo	Wax, yellow	6111407	1
Caelo	Wax, yellow	13335518	2
Caelo	Wax, yellow	14189905	1
Caelo	Wax, yellow	7091704	1
Caelo	Wax, yellow	13335503	6
Caelo	Wax, yellow	17011610	2
Caelo	Wax, yellow	4031403	1
Caelo	Wax, yellow	10121302	1
Caelo	Wax, yellow	15292022	1
Caelo	Wax, yellow	14189901	1
Caelo	Wax, yellow	11017318	3
Caelo	Wax, yellow	16034206	1
Caelo	Wax, yellow	11017314	1
Caelo	Wax, yellow	13031220	4
Caelo	Wax, yellow	13378630	2
Caelo	Wax, yellow	18299208	1
Caelo	Wax, yellow	15022907	1
Caelo	Wax, yellow	12367921	1
Caelo	Wax, yellow	12091404	1
Caelo	Wax, yellow	16038206	20
Caelo	Wax, yellow	15292005	4
Caelo	Wax, yellow	16038203	11
Caelo	Wax, yellow	15022902	6
Caelo	Wax, yellow	15023004	9
Caelo	Wax, yellow	12330116	1
Caelo	Wax, yellow	12367915	9
Caelo	Wax, yellow	18299201	1
Caelo	Wax, yellow	18113509	2
Caelo	Wax, yellow	16038207	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wax, yellow	16038205	2
Caelo	Wax, yellow	13335505	1
Caelo	Wax, yellow	17207404	2
Caelo	Wax, yellow	30011407	1
Caelo	Wax, yellow	17011403	1
Caelo	Wax, yellow	13335502	11
Caelo	Wax, yellow	17011612	20
Caelo	Wax, yellow	19031810	1
Caelo	Wax, yellow	17207405	2
Caelo	Wax, yellow	18049001	12
Caelo	Wax, yellow	18049003	3
Caelo	Wax, yellow	13031231	6
Caelo	Wax, yellow	1335503	1
Caelo	Wax, yellow	18049002	2
Caelo	Wax, yellow	18113501	19
Caelo	Wax, yellow	18113502	4
Caelo	Wax, yellow	18048906	2
Caelo	Wax, yellow	133335507	1
Caelo	Wax, yellow	18113506	2
Caelo	Wax, yellow	13378603	1
Caelo	Wax, yellow	15023006	3
Caelo	Wheat starch	11350907	1
Caelo	Wheat starch	1114611	1
Caelo	Wheat starch	18120004	1
Caelo	Wheat starch	16162905	1
Caelo	Wheat starch	15124707	1
Caelo	Wheat starch	1512477	1
Caelo	Wheat starch	15124710	1
Caelo	Wheat starch	15421507	2
Caelo	White deadnettle flower	11371613	2
Caelo	White deadnettle flower	15299101	1
Caelo	White deadnettle flower	12152402	2
Caelo	White deadnettle flower	12275802	1
Caelo	White deadnettle flower	11371610	5
Caelo	White deadnettle flower	52470437	1
Caelo	White deadnettle flower	2101301	1
Caelo	White deadnettle flower	11371611	2
Caelo	White deadnettle flower	10081206	1
Caelo	White deadnettle flower	11371609	3
Caelo	White deadnettle herb...	11371702	1
Caelo	White deadnettle herb...	11371703	1
Caelo	White deadnettle herb...	11371709	2
Caelo	White deadnettle herb...	11371711	2
Caelo	White deadnettle herb...	11371705	1
Caelo	White deadnettle herb...	11371712	1
Caelo	White deadnettle herb...	11371706	3
Caelo	White mustard seed, p...	13224804	3
Caelo	White mustard seed, p...	11114409	1
Caelo	White mustard seed, p...	13224801	1
Caelo	White mustard seed, p...	12375704	1
Caelo	White mustard seed, p...	12375706	2
Caelo	White mustard seed, p...	12375705	1
Caelo	White mustard seed, p...	12375705	4
Caelo	White mustard seed, p...	12375701	1
Caelo	White mustard seed, p...	12375702	1
Caelo	White mustard seed, p...	7254a110823	1
Caelo	Wild apple, sweet	11229805	1
Caelo	Wild pansy, cut	14061311	1
Caelo	Wild pansy, cut	12189203	5
Caelo	Wild pansy, cut	12189203	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Wild pansy, cut	11216103	2
Caelo	Wild pansy, cut	12382806	1
Caelo	Wild pansy, cut	13091205	1
Caelo	Wild pansy, cut	10291009	3
Caelo	Wild pansy, cut	17101205	1
Caelo	Wild pansy, cut	12382802	1
Caelo	Wild pansy, cut	1289203	1
Caelo	Wild pansy, cut	12189204	10
Caelo	Wild pansy, cut	12189204	4
Caelo	Wild pansy, cut	10291007	1
Caelo	Wild pansy, cut	12189204-3458312120	1
Caelo	Wild pansy, cut	15021301	1
Caelo	Wild pansy, cut	5021304	1
Caelo	Wild pansy, cut	12189204-3463312101	1
Caelo	Wild pansy, cut	7101303	1
Caelo	Wild thyme, cut	4501A120813	2
Caelo	Wild thyme, cut	11374805	5
Caelo	Wild thyme, cut	11374803	1
Caelo	Wild thyme, cut	11031203	1
Caelo	Wild thyme, cut	13086805	3
Caelo	Wild thyme, cut	11374805i.0.	1
Caelo	Willow bark, cut	13180705	1
Caelo	Willow bark, cut	21011302	1
Caelo	Willow bark, cut	13180705-3530411072	1
Caelo	Willow bark, cut	12185201	5
Caelo	Willow bark, cut	12185206	5
Caelo	Willow bark, cut	13180701	1
Caelo	Willow bark, cut	13180702	4
Caelo	Willow bark, cut	11073201	1
Caelo	Willow bark, cut	9071301	1
Caelo	Willow bark, cut	12185208	12
Caelo	Willow bark, cut	11073208	1
Caelo	Willow bark, cut	2185206	1
Caelo	Willow bark, cut	12185204	4
Caelo	Willow bark, cut	11073210	3
Caelo	Willow bark, cut	12185207	3
Caelo	Willow herb, cut	13104702	1
Caelo	Willow herb, cut	13200203	1
Caelo	Willow herb, cut	11347311	1
Caelo	Willow herb, cut	11347316	1
Caelo	Willow herb, cut	11347318	2
Caelo	Willowherb, small flo...	29121201	1
Caelo	Willowherb, small flo...	21061307	1
Caelo	Willowherb, small flo...	412041	1
Caelo	Willowherb, small flo...	2081305	1
Caelo	Willowherb, small flo...	10326108	14
Caelo	Willowherb, small flo...	12041301	1
Caelo	Willowherb, small flo...	10326107	1
Caelo	Willowherb, small flo...	1031301	1
Caelo	Willowherb, small flo...	13200201	3
Caelo	Willowherb, small flo...	14267909	1
Caelo	Willowherb, small flo...	12117701	8
Caelo	Willowherb, small flo...	31101303	1
Caelo	Willowherb, small flo...	13200203	7
Caelo	Willowherb, small flo...	9378908	1
Caelo	Willowherb, small flo...	13355906	2
Caelo	Willowherb, small flo...	10326104	4
Caelo	Willowherb, small flo...	15101302	1
Caelo	Willowherb, small flo...	28031318	1
Caelo	Willowherb, small flo...	2412E-01641	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Willowherb, small flo...	12117701-3371319012	1
Caelo	Willowherb, small flo...	13071301	1
Caelo	Willowherb, small flo...	7071201	1
Caelo	Willowherb, small flo...	12064506	1
Caelo	Willowherb, small flo...	16081304	1
Caelo	Willowherb, small flo...	16081302	1
Caelo	Willowherb, small flo...	18031305	1
Caelo	Willowherb, small flo...	13200202-3481312102	1
Caelo	Willowherb, small flo...	15011411	1
Caelo	Wool alcohols	15304505	2
Caelo	Wool alcohols	23101408	1
Caelo	Wool alcohols	4/2015	1
Caelo	Wool alcohols	13226408	1
Caelo	Wool alcohols	14262501	1
Caelo	Wool alcohols	13071502	1
Caelo	Wool alcohols	15304502	1
Caelo	Wool alcohols	6111801	1
Caelo	Wool alcohols	17169701	3
Caelo	Wool alcohols	18092502	2
Caelo	Wool alcohols	12340505	1
Caelo	Wool alcohols	13083906	2
Caelo	Wool alcohols	13083912	1
Caelo	Wool alcohols	11011401	1
Caelo	Wormwood herb, cut	12355804	8
Caelo	Wormwood herb, cut	272957	1
Caelo	Wormwood herb, cut	7031301	1
Caelo	Wormwood herb, cut	12355803	27
Caelo	Wormwood herb, cut	12355805	4
Caelo	Wormwood herb, cut	12355801	20
Caelo	Wormwood herb, cut	12355802	3
Caelo	Wormwood herb, cut	11285810	6
Caelo	Wormwood herb, cut	13381905	1
Caelo	Wormwood herb, cut	13174402	1
Caelo	Wormwood herb, cut	13174401	3
Caelo	Wormwood herb, cut	6101203	1
Caelo	Wormwood herb, cut	7011401	1
Caelo	Wormwood herb, cut	11285805	2
Caelo	Wormwood herb, cut	11285802	2
Caelo	Wormwood herb, cut	12355807	10
Caelo	Wormwood herb, cut	15408304	1
Caelo	Wormwood herb, cut	13381902	3
Caelo	Wormwood herb, cut	11285803	1
Caelo	Wormwood herb, cut	277513	1
Caelo	Wormwood herb, cut	12355808	2
Caelo	Wu Wei Zi	G220H0401421	3
Caelo	Xylitol	13230305	1
Caelo	Xylitol	5112I-02787	1
Caelo	Xylitol	13230309	1
Caelo	Xylitol	14068206	1
Caelo	Xylitol	14061401	1
Caelo	Xylitol	18113304	1
Caelo	Xylitol	16317401	1
Caelo	Xylitol	15388104	2
Caelo	Xylitol	12240612	1
Caelo	Xylitol	3314M-02787	1
Caelo	Xylitol	13230315	1
Caelo	Xylitol	14G04-B04-298356	2
Caelo	Xylitol	13230302	1
Caelo	Xylitol	18045302	4
Caelo	Xylometazoline hydroc...	16074205	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Xylometazoline hydroc...	17345706	1
Caelo	Xylometazoline hydroc...	17073905	6
Caelo	Xylometazoline hydroc...	13171202	11
Caelo	Xylometazoline hydroc...	13171203	3
Caelo	Xylometazoline hydroc...	12287212	2
Caelo	Xylometazoline hydroc...	160742095416	1
Caelo	Xylometazoline hydroc...	12287213	1
Caelo	Xylometazoline hydroc...	15037912	1
Caelo	Xylometazoline hydroc...	13344901	1
Caelo	Xylometazoline hydroc...	15292811	1
Caelo	Xylometazoline hydroc...	16074201	3
Caelo	Xylometazoline hydroc...	13344908	1
Caelo	Xylometazoline hydroc...	2031502	1
Caelo	Xylometazoline hydroc...	160740209	1
Caelo	Xylometazoline hydroc...	13171205	1
Caelo	Xylometazoline hydroc...	23051820	1
Caelo	Xylometazoline hydroc...	16074216	2
Caelo	Xylometazoline hydroc...	16074209	5
Caelo	Xylometazoline hydroc...	181102181	1
Caelo	Xylometazoline hydroc...	16074202	1
Caelo	Xylometazoline hydroc...	9278404	1
Caelo	Xylometazoline hydroc...	16074217	3
Caelo	Xylometazoline hydroc...	17345707	2
Caelo	Xylometazoline hydroc...	153037907	1
Caelo	Xylometazoline hydroc...	17073902	1
Caelo	Xylometazoline hydroc...	15292803	3
Caelo	Xylometazoline hydroc...	15037910	2
Caelo	Xylometazoline hydroc...	10118821	2
Caelo	Xylometazoline hydroc...	24111202	1
Caelo	Xylometazoline hydroc...	1703902	1
Caelo	Xylometazoline hydroc...	17345703	4
Caelo	Xylometazoline hydroc...	15292804	4
Caelo	Xylometazoline hydroc...	17345708	1
Caelo	Xylometazoline hydroc...	12287201	2
Caelo	Xylometazoline hydroc...	9276404	1
Caelo	Xylometazoline hydroc...	12287203	9
Caelo	Xylometazoline hydroc...	12287210	8
Caelo	Xylometazoline hydroc...	15037907	6
Caelo	Xylometazoline hydroc...	17073909	2
Caelo	Xylometazoline hydroc...	12287204	2
Caelo	Xylometazoline hydroc...	12027906	2
Caelo	Xylometazoline hydroc...	15292812	2
Caelo	Xylose	13133410	1
Caelo	Xylose	13133404	2
Caelo	Xylose	13133414	3
Caelo	Xylose	5031501	1
Caelo	Xylose	11031501	1
Caelo	Xylose	16002302	2
Caelo	Xylose	16002309	1
Caelo	Xylose	16206405	1
Caelo	Xylose	17160905	1
Caelo	Xylose	17160902	1
Caelo	Xylose	12120101	1
Caelo	Yarrow flower, rubbed	13350507	1
Caelo	Yarrow flower, rubbed	70139427	1
Caelo	Yarrow flower, rubbed	12267703	2
Caelo	Yarrow flower, rubbed	12363005	5
Caelo	Yarrow flower, rubbed	12267702	1
Caelo	Yarrow flower, rubbed	12363002	6
Caelo	Yarrow flower, rubbed	12338306	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Yarrow herb, cut	13029104	5
Caelo	Yarrow herb, cut	9121304	1
Caelo	Yarrow herb, cut	12379607	1
Caelo	Yarrow herb, cut	8011421	1
Caelo	Yarrow herb, cut	27011404	1
Caelo	Yarrow herb, cut	12379602	9
Caelo	Yarrow herb, cut	11280511	2
Caelo	Yarrow herb, cut	11280517	1
Caelo	Yarrow herb, cut	13381806	2
Caelo	Yarrow herb, cut	130929105	1
Caelo	Yarrow herb, cut	13381810	1
Caelo	Yarrow herb, cut	13381805	1
Caelo	Yarrow herb, cut	11280512	9
Caelo	Yarrow herb, cut	14018503	2
Caelo	Yarrow herb, cut	1408510	1
Caelo	Yarrow herb, cut	15006713	1
Caelo	Yarrow herb, cut	15006720	1
Caelo	Yarrow herb, cut	15006708	1
Caelo	Yarrow herb, cut	12338305	7
Caelo	Yarrow herb, cut	12338303	10
Caelo	Yarrow herb, cut	5021302	1
Caelo	Yarrow herb, cut	12379606	9
Caelo	Yarrow herb, cut	13029107-3436311151	1
Caelo	Yarrow herb, cut	13029107-3418310182	1
Caelo	Yarrow herb, cut	H-220313-1	1
Caelo	Yarrow herb, cut	10101205	1
Caelo	Yarrow herb, cut	12379601	9
Caelo	Yarrow herb, cut	16041302	1
Caelo	Yarrow herb, cut	12041302	1
Caelo	Yarrow herb, cut	29071301	1
Caelo	Yarrow herb, cut	10081301	1
Caelo	Yarrow herb, cut	13029108	9
Caelo	Yarrow herb, cut	12338311	2
Caelo	Yarrow herb, cut	13029106	10
Caelo	Yarrow herb, cut	13029103	12
Caelo	Yarrow herb, cut	12338304	15
Caelo	Yarrow herb, cut	12338306-3358319071	1
Caelo	Yarrow herb, cut	9101311	1
Caelo	Yarrow herb, cut	12336306	1
Caelo	Yarrow herb, cut	11306301	1
Caelo	Yarrow herb, cut	16102013E	1
Caelo	Yarrow herb, cut	1309106	1
Caelo	Yarrow herb, cut	11280514	2
Caelo	Yarrow herb, cut	12379603	7
Caelo	Yarrow herb, cut	11218404	1
Caelo	Yarrow herb, cut	13029105	4
Caelo	Yellow chaste weed, cut	15029505	1
Caelo	Yellow chaste weed, cut	11293101	4
Caelo	Yellow chaste weed, cut	10244407	1
Caelo	Yellow chaste weed, cut	512U-01473	1
Caelo	Yellow chaste weed, cut	12170602	2
Caelo	Yellow chaste weed, cut	12170605	2
Caelo	Yellow chaste weed, cut	12170603	1
Caelo	Yellow chaste weed, cut	1024407	1
Caelo	Zedoary root, cut	25111304	1
Caelo	Zedoary root, cut	12125304	2
Caelo	Zedoary root, cut	13010104	1
Caelo	Zedoary root, cut	13010102	3
Caelo	Zedoary root, cut	11091603	1
Caelo	Zedoary root, cut	111091607	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caelo	Zedoary root, cut	12125301	3
Caelo	Zedoary root, powdered	13165504	4
Caelo	Zedoary root, powdered	15277101	1
Caelo	Zedoary root, powdered	15352606	1
Caelo	Zedoary root, powdered	12007504	1
Caelo	Zedoary root, powdered	13165501	2
Caelo	Zinc acetate dihydrate	15383306	3
Caelo	Zinc acetate dihydrate	18204604	1
Caelo	Zinc acetate dihydrate	12134803	1
Caelo	Zinc acetate dihydrate	15160801	1
Caelo	Zinc acetate dihydrate	15383308	1
Caelo	Zinc acetate dihydrate	12134801	2
Caelo	Zinc acetate dihydrate	12316405	1
Caelo	Zinc acetate dihydrate	17197601	5
Caelo	Zinc acetate dihydrate	12316403	2
Caelo	Zinc acetate dihydrate	17020404	2
Caelo	Zinc acetate dihydrate	12316404	2
Caelo	Zinc acetate dihydrate	17365807	1
Caelo	Zinc acetate dihydrate	17020403	2
Caelo	Zinc acetate dihydrate	15383304	1
Caelo	Zinc acetate dihydrate	12316401	2
Caelo	Zinc acetate dihydrate	15160804	1
Caelo	Zinc gluconate	17095103	1
Caelo	Zinc gluconate	16178801	1
Caelo	Zinc gluconate	17016502	2
Caelo	Zinc sulfate heptahyd...	12257816	1
Caelo	Zinc sulfate heptahyd...	11095601	1
Caelo	Zinc sulfate heptahyd...	12257818	1
Caelo	Zinc sulfate heptahyd...	4012Q-02798	1
Caelo	Zinc sulfate heptahyd...	12257802	1
Caelo	Zinc sulfate heptahyd...	12257809	2
Caelo	Zinc sulfate heptahyd...	81554199	1
Caelo	Zinc sulfate heptahyd...	11095621	1
Caeo	Glucose monohydrate	1049103	1
Caeo	Mometasone furoate, m...	12351201	1
Caeo/Noweda	Lactose monohydrate	1523204	1
Caeo/Noweda	Mometasone furoate, m...	15279003	1
Caeolo	Emulsifying cetostear...	17056807	1
Caeolo/Gehe	Magnesium hydrogen ci...	17173405	2
Caeolo/Gehe	Menthol	16096619	1
Caeolo/Krieger	Erythromycin	17088004	9
Caeolo/Noweda	Glucose monohydrate	15413301	2
Caeolo/Phönix	Triclosan	18081416	8
Caesar	Mannitol	1508012-1	1
Caesar & Lorenz	Yarrow herb, cut	12379605	3
Caesar & Lorenz GmbH	Ciclopirox olamine	171188007	1
Caesar & Loretz	Citric acid	15087904	3
Caesar & Loretz	Hydrocortisone	17256603	8
Caesar & Loretz	Hydrocortisone acetate	18142815	4
Caesar & Loretz	Stearic acid 50	14359310	1
Caesar & Loretz GmbH	Acetylsalicylic acid	19116604	1
Caesar & Loretz GmbH	Acetylsalicylic acid	15334914	4
Caesar & Loretz GmbH	Allantoin	12374403	10
Caesar & Loretz GmbH	Allantoin	12374416	1
Caesar & Loretz GmbH	Betamethasone valerate	18199918	1
Caesar & Loretz GmbH	Bismuth gallath, basic	17259510	2
Caesar & Loretz GmbH	Bismuth gallath, basic	17259513	1
Caesar & Loretz GmbH	Cetostearyl alcohol	17056806	1
Caesar & Loretz GmbH	Cetostearyl alcohol	18248502	1
Caesar & Loretz GmbH	Cetyl palmitate 15	10068402	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz GmbH	Cetyl palmitate 15	18149912	2
Caesar & Loretz GmbH	Chloral hydrate	17314317	1
Caesar & Loretz GmbH	Chloral hydrate	17314319	1
Caesar & Loretz GmbH	Chlortetracycline hyd...	12C16-N01	2
Caesar & Loretz GmbH	Ciclopirox olamine	17307410	1
Caesar & Loretz GmbH	Citric acid	18129207	1
Caesar & Loretz GmbH	Clobetasol propionate	18082709	7
Caesar & Loretz GmbH	Clobetasol propionate	18082710	3
Caesar & Loretz GmbH	Clotrimazole	18147216	2
Caesar & Loretz GmbH	Clotrimazole	19016105	2
Caesar & Loretz GmbH	Cocoa butter	12270210	3
Caesar & Loretz GmbH	Dexamethasone	15054622	26
Caesar & Loretz GmbH	Dithranol	19086804	1
Caesar & Loretz GmbH	Dithranol	19055001	1
Caesar & Loretz GmbH	Emulsifying cetostear...	18248502	1
Caesar & Loretz GmbH	Emulsifying cetostear...	16119604	11
Caesar & Loretz GmbH	Emulsifying cetostear...	18291403	1
Caesar & Loretz GmbH	Emulsifying cetostear...	17132001	6
Caesar & Loretz GmbH	Ephedrine hydrochloride	14378503	1
Caesar & Loretz GmbH	Erythromycin	18096205	1
Caesar & Loretz GmbH	Estradiol valerate	18067415	1
Caesar & Loretz GmbH	Estriol	19035609	1
Caesar & Loretz GmbH	Fructose	18291602	1
Caesar & Loretz GmbH	Fructose	16283910	6
Caesar & Loretz GmbH	Gentamicin sulfate	17358607	1
Caesar & Loretz GmbH	Glucose	18012926	1
Caesar & Loretz GmbH	Glucose monohydrate	15180013	3
Caesar & Loretz GmbH	Hawthorn flower	11068206	3
Caesar & Loretz GmbH	Hydrocortisone	16004612	11
Caesar & Loretz GmbH	Hydrocortisone	18091526	2
Caesar & Loretz GmbH	Lactose monohydrate	18198006	1
Caesar & Loretz GmbH	Lactose monohydrate	14140311	1
Caesar & Loretz GmbH	Lidocaine	18044414	1
Caesar & Loretz GmbH	Lidocaine hydrochlori...	13343002	4
Caesar & Loretz GmbH	Lidocaine hydrochlori...	15215201	4
Caesar & Loretz GmbH	Macrogols 4000	19009313	1
Caesar & Loretz GmbH	Macrogols 4000	14180520	1
Caesar & Loretz GmbH	Macrogols 4000	18315307	3
Caesar & Loretz GmbH	Macrogols 4000	18244806	1
Caesar & Loretz GmbH	Magnesium chloride he...	13064223	1
Caesar & Loretz GmbH	Mannitol	14323119	1
Caesar & Loretz GmbH	Methyl-4-hydroxybenzo...	18176910	1
Caesar & Loretz GmbH	Methylcellulose	18331803	2
Caesar & Loretz GmbH	Methylcellulose	12044922	3
Caesar & Loretz GmbH	Metronidazole	19010201	2
Caesar & Loretz GmbH	Mometasone furoate, m...	18334101	2
Caesar & Loretz GmbH	Octenidine dihydrochl...	15006210	1
Caesar & Loretz GmbH	Octenidine dihydrochl...	181975002	1
Caesar & Loretz GmbH	Octenidine dihydrochl...	18223423	1
Caesar & Loretz GmbH	Paracetamol	18083008	1
Caesar & Loretz GmbH	Paracetamol	18083006	1
Caesar & Loretz GmbH	Peppermint leaf, cut	11294411	1
Caesar & Loretz GmbH	Poloxamers 407	15183102	3
Caesar & Loretz GmbH	Potassium sorbate	11261303	1
Caesar & Loretz GmbH	Povidone 25, soluble	18300002	1
Caesar & Loretz GmbH	Prednicarbate, micron...	19051602	1
Caesar & Loretz GmbH	Prednisolone acetate	18067314	2
Caesar & Loretz GmbH	Procaine hydrochloride	14164413	2
Caesar & Loretz GmbH	Progesterone, microni...	17348409	1
Caesar & Loretz GmbH	Saccharin sodium	17329311	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz GmbH	Salicylic acid	16312012	2
Caesar & Loretz GmbH	Sodium citrate	14369406	6
Caesar & Loretz GmbH	Sodium citrate	13089334	1
Caesar & Loretz GmbH	Sodium cyclamate	12105622	1
Caesar & Loretz GmbH	Softisan [®] 601	18197002	1
Caesar & Loretz GmbH	Soy lecithin (instant...	14281603	1
Caesar & Loretz GmbH	Spironolactone	12049129	2
Caesar & Loretz GmbH	Testosterone propionate	17325717	2
Caesar & Loretz GmbH	Tetracaine base	18124906	1
Caesar & Loretz GmbH	Thiomersal	12275218	1
Caesar & Loretz GmbH	Thymol	17340716	1
Caesar & Loretz GmbH	Tretinoin	11225515	1
Caesar & Loretz GmbH	Tretinoin	12320211	1
Caesar & Loretz GmbH	Triamcinolone acetonide	17057906	10
Caesar & Loretz GmbH	Urea	15087105	12
Caesar & Loretz GmbH	Urea	1306421	1
Caesar & Loretz GmbH	Wax, bleached	18239104	2
Caesar & Loretz GmbH	Wax, yellow	13031221	2
Caesar & Loretz GmbH ...	Hydroxyethylcellulose...	12274702	2
Caesar & Loretz GmbH ...	Bitter fennel	11016331	1
Caesar & Loretz GmbH ...	Potassium hydroxide	12057902	2
Caesar & Loretz GmbH ...	Pyridoxine hydrochlor...	11197009	6
Caesar & Loretz GmbH ...	Sage leaf, cut	11056121	3
Caesar & Loretz GmbH ...	Salicylic acid	11311904	1
Caesar & Loretz GmbH ...	Vanillin	10066830	2
Caesar & Loretz GmbH ...	Saccharin sodium	13205418	3
Caesar & Loretz GmbH,...	Aluminium chloride he...	12323702	14
Caesar & Loretz GmbH,...	Aluminium chloride he...	12047508	4
Caesar & Loretz GmbH,...	Betamethasone valerate	10296707	2
Caesar & Loretz GmbH,...	Birch leaf, cut	11050102	3
Caesar & Loretz GmbH,...	Bitter fennel	11016319	1
Caesar & Loretz GmbH,...	Camphor	11032035	4
Caesar & Loretz GmbH,...	Caraway fruit	11003213	2
Caesar & Loretz GmbH,...	Caraway fruit	9267824	2
Caesar & Loretz GmbH,...	Cetostearyl alcohol	11302101	2
Caesar & Loretz GmbH,...	Cyanocobalamin	9324102	2
Caesar & Loretz GmbH,...	Dexamethasone	10369922	1
Caesar & Loretz GmbH,...	Doxylamine hydrogen s...	13A04-N02	7
Caesar & Loretz GmbH,...	Ephedrine hydrochloride	11126411	1
Caesar & Loretz GmbH,...	Equisetum stem, cut	10375510	2
Caesar & Loretz GmbH,...	Equisetum stem, cut	12288801	16
Caesar & Loretz GmbH,...	Erythromycin	10217208	1
Caesar & Loretz GmbH,...	Gentian root, cut	12244603	2
Caesar & Loretz GmbH,...	Liquorice root, peeled	11034014	1
Caesar & Loretz GmbH,...	Lungwort herb, cut	13029304	5
Caesar & Loretz GmbH,...	Mallow flower, cut	13267301	5
Caesar & Loretz GmbH,...	Mallow flower, cut	12187402	7
Caesar & Loretz GmbH,...	Maté leaf, cut	13118002	10
Caesar & Loretz GmbH,...	Menthol	11202910	6
Caesar & Loretz GmbH,...	Menthol	12020402	6
Caesar & Loretz GmbH,...	Miconazole nitrate	11324904	1
Caesar & Loretz GmbH,...	Neomycin sulfate	11/19-N14	5
Caesar & Loretz GmbH,...	Nettle leaf, cut	10291423	6
Caesar & Loretz GmbH,...	Pyridoxine hydrochlor...	12214119	4
Caesar & Loretz GmbH,...	Pyridoxine hydrochlor...	12214115	7
Caesar & Loretz GmbH,...	Pyridoxine hydrochlor...	12214102	7
Caesar & Loretz GmbH,...	Pyridoxine hydrochlor...	12214110	2
Caesar & Loretz GmbH,...	Raspberry leaf, cut	10013408	1
Caesar & Loretz GmbH,...	Ribwort plantain, cut	12362903	29
Caesar & Loretz GmbH,...	Ribwort plantain, cut	13113004	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar & Loretz GmbH,...	Sage leaf, cut	12370016	8
Caesar & Loretz GmbH,...	Sage leaf, cut	11347502	5
Caesar & Loretz GmbH,...	Sage leaf, cut	12370010	16
Caesar & Loretz GmbH,...	Sage leaf, cut	12370002	4
Caesar & Loretz GmbH,...	Selenium yeast	11318111	1
Caesar & Loretz GmbH,...	Senna leaf, cut	11277314	5
Caesar & Loretz GmbH,...	Senna leaf, cut	11277308	3
Caesar & Loretz GmbH,...	Senna leaf, cut	13136708	7
Caesar & Loretz GmbH,...	Senna leaf, cut	11277303	3
Caesar & Loretz GmbH,...	Senna leaf, cut	12273103	11
Caesar & Loretz GmbH,...	Senna pods, tinnevelly	12267603	1
Caesar & Loretz GmbH,...	Senna pods, tinnevelly	12267601	2
Caesar & Loretz GmbH,...	Senna pods, tinnevelly	11003014	1
Caesar & Loretz GmbH,...	Senna pods, tinnevelly	11003012	1
Caesar & Loretz GmbH,...	Sodium edetate	12181711	11
Caesar & Loretz GmbH,...	Sodium hydrogen carbo...	13060102	1
Caesar & Loretz GmbH,...	Soy lecithin (instant...	12139013	1
Caesar & Loretz GmbH,...	Tall golden rod, cut	11207014	3
Caesar & Loretz GmbH,...	Trometamol	11160714	3
Caesar & Loretz GmbH,...	Urea	11150736	5
Caesar & Loretz GmbH,...	Urea	12111828	22
Caesar & Loretz GmbH,...	Urea	12111804	5
Caesar & Loretz GmbH,...	Vanillin	12079412	4
Caesar & Loretz GmbH,...	Vanillin	12079401	1
Caesar & Loretz GmbH,...	Wax, yellow	11285008	1
Caesar & Loretz GmbH,...	Wild thyme, cut	11374802	1
Caesar & Loretz GmbH,...	Yarrow herb, cut	12338306	3
Caesar & Loretz GmbH,...	Prednisolone acetate	17185601	18
Caesar & Loretz GmbH/...	Macrogols 4000	19009302	1
Caesar & Loretz GmbH/...	Sodium citrate	18208315	2
Caesar & Loretz GmbH/...	Trometamol	18197108	1
Caesar & Loretz GmbH/...	Ethacridine lactate m...	17332316	2
Caesar & Loretz GmbH/...	Cocoa butter	18254402	1
Caesar & Loretz GmbH/Ph	Estradiol valerate	18067412	2
Caesar & Loretz GmbH/...	Lidocaine	18044410	1
Caesar & Loretz GmbH/...	Estradiol benzoate	18144404	5
Caesar & Loretz/ Caes...	Glucose monohydrate	18030004	6
Caesar & Loretz/ Noweda	Heavy alkaline bismut...	12145701	1
Caesar & Loretz/ Noweda	Hyssop herb, cut	11374301	1
Caesar & Loretz/ Noweda	Ivy leaf, cut	12225606	5
Caesar & Loretz/Noweda	Aniseed	153986002	1
Caesar & Loretz/Noweda	Caraway fruit	15056601	1
Caesar & Loretz GmbH	Emulsifying cetostear...	12065703	1
Caesar & Loretz GmbH	Soy lecithin (instant...	13279813	1
Caesar § Lorenz GmbH	Estriol	12298712	1
Caesar § Lorenz GmbH	Thymol	13288806	1
Caesar Loretz	Camphor	18327903	1
Caesar u Loretz	Cetyl palmitate 15	14213710	2
Caesar u Loretz	Clotrimazole	16234501	9
Caesar u Loretz	Propranolol hydrochlo...	13379807	5
Caesar u Loretz	Spironolactone	12049119	2
Caesar u. Loretz/Noweda	Erythromycin	16079209	2
Caesar u. Loretz/Noweda	Prednisolone, microni...	16134301	5
Caesar u.Loretz	Iodoform	16012209	2
Caesar&Lorenz	Clobetasol propionate	1610370	1
Caesar&Lorenz	Progesterone, microni...	16103605	3
Caesar&Lorenz/Phoenix	Alfatradiol	16087502	8
Caesar&Loretz	Camphor	14310719	1
Caesar&Loretz	Captopril	17210201	2
Caesar&Loretz GmbH D...	Miconazole nitrate	12243808	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar&Loretz GmbH	Erythromycin	15042807	1
Caesar&Loretz GmbH	Potassium sorbate	17170014	18
Caesar&Loretz GmbH	Prednisolone acetate	15396408	1
Caesar&Loretz GmbH /A...	Lidocaine hydrochlori...	15215303	3
Caesar&Loretz GmbH D-...	Senna leaf, cut	12115504	4
Caesar&Loretz GmbH D-...	Urea	12111808	14
caesar&Loretz GmbH,Ph...	Mometasone furoate, m...	16028604	6
Caesar&Loretz GmbH/ A...	Oxalic acid	13353920	1
Caesar&Loretz GmbH/An...	Lidocaine hydrochlori...	13118803	7
Caesar&Loretz GmbH/An...	Urea	13064217	2
Caesar&Loretz,AHC	Nystatin	14196509	2
Caesar&Loretz,AHC,08...	Mometasone furoate, m...	16028601, f (NRF) 1,014	1
Caesar&Loretz,Anzag	Nystatin	12158008	7
Caesar&Loretz,Anzag	Sodium monohydrogenph...	12256606	3
Caesar&Loretz,Anzag	Urea	11150738	2
Caesar&Loretz,Phoenix	Erythromycin	13175602	3
Caesar&Loretz,Phoenix	Frangula bark, cut	14163104	1
Caesar&Loretz,Phoenix	Hydrocortisone	15181004	9
Caesar&Loretz,Phoenix	Ladys bedstraw herb,...	14019206	1
Caesar&Loretz,Phoenix	Ladys bedstraw herb,...	14019215	1
Caesar&Loretz,Phoenix	Marjoram herb, rubbed	14181005	1
Caesar&Loretz,Phoenix	Methyl-4-hydroxybenzo...	13268821	1
Caesar&Loretz,Phoenix	Metoclopramide hydroc...	14193101	2
Caesar&Loretz,Phoenix	Metronidazole	12110410	1
Caesar&Loretz,Phoenix	Neomycin sulfate	12C13-N04	5
Caesar&Loretz,Phoenix	Nystatin	1218008	1
Caesar&Loretz,Phoenix	Nystatin	1111057-01	8
Caesar&Loretz,Phoenix	Nystatin	13173408	2
Caesar&Loretz,Phoenix	Nystatin	14363512	1
Caesar&Loretz,Phoenix	Salicylic acid	11197104	2
Caesar&Loretz,Phoenix	Salicylic acid	14076212	2
Caesar&Loretz,Phoenix	Salicylic acid	9220710	1
Caesar&Loretz,Phoenix	Spironolactone	11042603	1
Caesar&Loretz,Phoenix	Spironolactone	11042613	1
Caesar&Loretz,Phoenix	Theophylline	12166315	2
Caesar&Loretz,Phoenix	Triamcinolone acetone...	12301502	1
Caesar&Loretz,Phoenix	Triamcinolone acetone...	11267705	2
Caesar&Loretz,Phoenix	Triamcinolone acetone...	12173703	1
Caesar&Loretz,Phoenix	Urea	11150713	2
Caesar&Loretz,Phoenix...	Prednisolone acetate	15398302, f (NRF) 0,98	1
caesar&Loretz,phönix	Betamethasone valerate	14188606	2
caesar&Loretz,phönix	Erythromycin	15042603	7
caesar&Loretz,Phönix	Iceland moss, cut	14304404	3
Caesar&Loretz,Phönix	Triclosan	15203606	10
Caesar&Loretz/Caesar&...	Urea	15154104	83
Caesar&Loretz/Noweda	Neomycin sulfate	15298904	1
Caesar/ Noweda	Black mustard seed, p...	14293711	3
Caesar/Lorenz	Clobetasol propionate	12342305	4
Caesar/Lorenz	Emulsifying cetosteary...	13233005	11
Caesar/Lorenz	Triclosan	13247904	14
Caesar/Noweda	Aloes, cape, cut	2013E-01089	1
Caesar/Noweda	Aloes, cape, cut	12338211	1
Caesar/Noweda	Aluminium potassium s...	12184411	2
Caesar/Noweda	American liverwort, cut	11254608	1
Caesar/Noweda	Bogbean leaf, cut	1510U-01517	1
Caesar/Noweda	Camphor	1032	1
Caesar/Noweda	Camphor	13195328	1
Caesar/Noweda	Cetostearyl alcohol	15272705	1
Caesar/Noweda	Cetostearyl alcohol	2911E-01077	1
Caesar/Noweda	Cetostearyl alcohol	2911E-01077	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Caesar/Noweda	Cetyl palmitate 15	13295909	1
Caesar/Noweda	Cetyl palmitate 15	12228801	1
Caesar/Noweda	Cetyl palmitate 15	12124811	1
Caesar/Noweda	Coriander	2010M-01539	1
Caesar/Noweda	Couch grass rhizome, ...	12186813	4
Caesar/Noweda	Diphenhydramine hydro...	12007714	1
Caesar/Noweda	Elder flower, rubbed	4910A-01471	1
Caesar/Noweda	Emulsifying cetostear...	11343405	1
Caesar/Noweda	Emulsifying cetostear...	12065707	1
Caesar/Noweda	Glycerol monostearate...	13221311	2
Caesar/Noweda	Heather, cut	4911I-03224	1
Caesar/Noweda	Ispaghula seed	12276516	1
Caesar/Noweda	Ispaghula seed	12276503	1
Caesar/Noweda	Ispaghula seed	12276509	3
Caesar/Noweda	Lycopodium herb, cut	1255008	1
Caesar/Noweda	Menthol	15076111	2
Caesar/Noweda	Menthol	12020403	1
Caesar/Noweda	Mullein flower, cut	2810E-01475	1
Caesar/Noweda	Propyl-4-hydroxybenzo...	3710M-01931	2
Caesar/Noweda	Psyllium seed	11185817	3
Caesar/Noweda	Psyllium seed	12276501	3
Caesar/Noweda	Restharrow root, cut	71015/3	2
Caesar/Noweda	Sodium citrate	12045008	1
Caesar/Noweda	Theriac, powdered	12289309	1
Caesar/Noweda	Urea	13064208	1
Caesar/Noweda	Wax, yellow	12124811	3
Caesar/Phoenix	Chloral hydrate	17250601	6
Caleo	Dexamethasone	15054602	8
Caleo/Phoenix	Tannin	14K17-B08-306993	1
Caleo/Phönix	Betamethasone valerate	1522609	1
Calo/AHD	Camphor	15086803	1
Calo/AHD	Erythromycin	15042605	7
Calo/AHD	Hyaluronic acid sodiu...	14364704	1
Camphor & Allied Prod...	Camphor	16209302	1
Canea Pharma GmbH	Nystatin	18072121	4
canea, gehe	Ascorbic acid	VC1242	2
canea, gehe	Potassium sorbate	13065215	1
carlo	Alfatradiol	10319903	1
cealo	Polymyxin B sulfate	16c01-b04-319432	1
Cealo	Betamethasone valerate	17078716	2
Cealo	Calcium citrate tetra...	12310401	2
Cealo	Chlorhexidine diacetate	13177002	12
Cealo	Ciclopirox olamine	14147403	2
Cealo	Clotrimazole	17111612	2
Cealo	Glucose monohydrate	13157017	1
Cealo	Ispaghula seed	12276502	4
Cealo	Lactose monohydrate	13104901	9
Cealo	Magnesium chloride he...	112354615	1
Cealo	Magnesium chloride he...	11354615	1
Cealo	Magnesium hydrogen ci...	13233907	4
Cealo	Magnesium hydrogen ci...	12368014	3
Cealo	Melatonin	12321422	1
Cealo	Pheniramine maleate	12352608	1
Cealo	Potassium sorbate	1710013	1
Cealo	Ribwort plantain, cut	12362907	1
Cealo	Senna leaf, cut	12273116	15
Cealo	Triclosan	71215/1Gie	1
Cealo	Triclosan	17081723	6
Cealo/ Noweda	Camphor	13108116	1
Cealo/ Noweda	Saccharin sodium	17329407	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Cealo/ Noweda	Wool alcohols	13380601	2
Cealo/AHD	Clindamycin hydrochlo...	16D05-B08-325088	4
Cealo/AHD	Paraffin, hard	15295411	2
Cealo/Alliance Health...	Cetostearyl alcohol	18291403	1
Cealo/Alliance Health...	Nystatin	18072118	1
Cealo/Anzag	Mometasone furoate, m...	18026012	2
Cealo/Gehe	Ambroxol hydrochloride	171158204	1
Cealo/Gehe	Betamethasone valerate	1739002	2
Cealo/Gehe	Macrogols 4000	16023907	1
Cealo/Gehe	Melatonin	10296624	1
Cealo/Gehe	Tormentil	14140006	3
Cealo/Noweda	Prednisolone, microni...	1526108	1
Cealo/Sanacorp	Captopril	12339617	2
Cealo/Sanacorp	Corn starch	17287302	2
Cealo/Sanacorp	Octenidine dihydrochl...	15006207	6
Cealo/Sanacorp	Octenidine dihydrochl...	16325305	10
Cealo/Sanacorp	Prednicarbate, micron...	15546306	1
Cesar & Loretz GmbH	Chlorhexidine diacetate	80104278	1
Celo	Agrimony, cut	11374003	1
Celo	Prednicarbate, micron...	1687408	1
Celo/Kehr	Tretinoin	16217122	3
cesar&lorentz/Phönix	Triamcinolone acetone...	13124801	5
Changzhou Yabang Phar...	Propranolol hydrochlo...	15F16-B02	2
Chem & Pharma	Ambroxol hydrochloride	AXL0361201310	1
Chiesi/ Now	Clobetasol propionate	13269503	3
Chiesi/Phönix	Black cumin seed	14359906	2
China Medica	Hard fat in pastilles	3244303	1
China Medica	Procaine hydrochloride	13109911	1
China Medica	Raspberry leaf, cut	13086703	1
Chinosol/Anzag	Ethacridine lactate m...	2112A-02196	1
Chinosol/Anzag	Ethacridine lactate m...	914E-02196	1
Chiracon GmbH	Alfatradiol	AP35-17-II-85P3624	1
Clariant - Mindel-Food	Bentonite	1301	1
Concord Biotech Ltd./...	Ciclosporin	14G22-B02-296891	1
Cranach-Apotheke	Sodium edetate	11a31-no6	1
Creapure	Creatine	515941	1
Crystal	Finasteride (form III)	141100801	1
Crystal Pharma	Clobetasol propionate	C0054/1	1
crystal pharma/gehe	Budesonide, micronized	1305004-01	2
Dan Pharma	Urea	D08118	1
Defektur	Bitter fennel	7/13	1
Defektur	Bitter fennel	57/13	1
Defektur	Bitter fennel	52/13	1
Defektur	Matricaria flower	10/13	1
Defektur	Matricaria flower	59/13	1
Defektur	Matricaria flower	29/13	1
Defektur	Matricaria flower	56/12	1
Defektur/eigen	Mannitol	1604008-01	34
Degussa / falcento	Creatine	422031	1
Delta	Betamethasone dipropi...	18B06-B05-346219	3
Denk	Levocarnitine	A04217	1
Denk	Pyridoxal-5-phosphate	C041402802	1
Denk	Pyridoxal-5-phosphate	20171202	1
Dermapharm	8-quinolinol sulfate ...	120905	1
Dermapharm	8-quinolinol sulfate ...	100804	1
Dermapharm	8-quinolinol sulfate ...	140804	1
Dermapharm	Ethacridine lactate m...	130201	7
Dermapharm	Ethacridine lactate m...	130201	2
Dermapharm	Ethacridine lactate m...	111201	7
Dermapharm	Quinolinol sulfate po...	100804	12

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Dermapharm	Quinolinol sulfate po...	120905	42
Dermapharm	Quinolinol sulfate po...	140804	71
Dermapharm	Quinolinol sulfate po...	1470804	1
Dermopharm	Ethacridine lactate m...	141101	11
Dextro Energy	Glucose monohydrate	13171	1
Dextro Energy GmbH/Gehe	Glucose monohydrate	L4163	2
Dextro Energy GmbH/Gehe	Glucose monohydrate	L4274	1
Dextro Energy GmbH/Gehe	Glucose monohydrate	L4164	1
Dextro Energy GmbH/Sa...	Glucose monohydrate	L1299	2
Dextro Energy GmbH	Glucose monohydrate	L4273	1
Dextro Energy GmbH	Glucose monohydrate	L4036	1
Dextro GmbH/Noweda	Glucose monohydrate	L5301	1
Dextro/Krieger/15.05.13	Glucose monohydrate	13092	1
DHU	Sucrose	9050610	1
DHU	Sucrose	9010710	1
DHU	Sucrose	9010709	1
DHU	Sucrose	9020110	1
DHU	Sucrose	7020407	1
DHU	Sucrose	9020510	1
Diamant	Sucrose	L24223	1
DIV GH	Macrogols 4000	13/12	1
DOW Chemical Company ...	Hypromellose 4000	14G02-B11-296850	7
Dr. August Wolff	Urea	13E17N02	1
Dr. Behr	Dimercaptosuccinic acid	WE14020002	1
Dr. Behr	Ubidecarenone	WE18020018	1
Dr. Paul Lohmann	Magnesium hydrogen ci...	1078434	1
Dr. Paul Lohmann	Magnesium hydrogen ci...	1077866	2
Dr. Paul Lohmann	Magnesium hydrogen ci...	1079474	1
Dr.Hollborn	Vanillin	6082015	1
Dr.Lohmann	Calcium citrate tetra...	3042527	1
Dr.Wolff	Urea	1234567	1
DSM Nutritional Prod....	Nicotinamide	1106044-01	1
ebe a	Isosorbide dinitrate ...	14D09-B03-294281	7
Ebert 24.08.16	Fusidic acid	1510-B02_319452	1
Ebert & Jacobi 12....	Citric acid	1506214	1
Ebert & Jacobi 12....	Xylose	16002301	3
Ebert KKH 10.10.16	Chloral hydrate	15448308	1
ebert u.jacobi	Marian thistle fruit	287527	1
Ebert-Jacobi	Cetostearyl alcohol	16230603	3
Ebert-Jacobi	Cholesterol	16044408	3
Ebert-Jacobi	Cyproterone acetate	16K09-B04-329036	1
Ebert-Jacobi	Quinolinol sulfate po...	16221804	2
Ebert-Jacobi	Tetracaine hydrochlor...	15K04-B03-328163	1
Edelweiss	Lactose monohydrate	253292176-1124	1
Eigen/Defektur	Triamcinolone acetone	U0544/113530	1
Eigenherstellung	Betamethasone valerate	9070313	1
Eigenherstellung	Glucose monohydrate	30912HML01	1
Eigenherstellung	Glucose monohydrate	110913HJK02	1
Eigenherstellung	Glucose monohydrate	40615HJK01	1
Eigenherstellung	Glucose monohydrate	160915HML01	1
Eigenherstellung	Glucose monohydrate	110912HML01	1
Eigenherstellung	Glucose monohydrate	120912HML02	1
Eigenherstellung	Glucose monohydrate	180213HJK01	1
Eigenherstellung	Lactose monohydrate	40113HJK01	1
Eigenherstellung	Lactose monohydrate	70113HML01	1
Eigenherstellung	Lactose monohydrate	110913HJK01	1
Eigenherstellung	Lactose monohydrate	60814HSA01	1
Eigenherstellung	Lactose monohydrate	40815HMS01	1
Eigenherstellung	Methyl-4-hydroxybenzo...	100912PML01	1
Eigenherstellung	Xylose	80414HJK01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
el compra import	Taurine	1609025	1
elcompra	Taurine	xp17013244	1
eu	Gentamicin sulfate	1610038-01	1
Eu Rho	Citric acid	L1202001-01	6
Eu Rho	Magnesium hydrogen ci...	110600401	1
Eu Rho	Magnesium hydrogen ci...	1202034-01	1
Eu Rho	Sodium hydrogen carbo...	1203030-01	8
Eurho	Chloramphenicol	1704087-03	1
Eurho	Fluocinolonacetomid, ...	1408022-04	1
Eurho	Menthol	170502201	4
Eurho	Triamcinolone acetamide	1805001-02	9
EuRho	alpha-tocopheryl acet...	1810024-01	1
EuRho	Aluminium chloride he...	1106017-01	1
EuRho	Aluminium chloride he...	2512I01095	1
EuRho	Aluminium chloride he...	d024a	1
EuRho	Ambroxol hydrochloride	1602015-02	4
EuRho	Betamethasone valerate	181008-01	1
EuRho	Betamethasone valerate	180102401	2
EuRho	Caffeine	17030006-02	2
EuRho	Calcium citrate tetra...	1502012-02	1
EuRho	Calcium gluconate	12090006-01	1
EuRho	Camphor	2111E-01219	1
EuRho	Capsaicin	1703023-02	2
EuRho	Cellulose, microcryst...	1201004-02	3
EuRho	Chloramphenicol	216A-01259	1
EuRho	Citric acid	L1503027-02	5
EuRho	Clioquinol	1309009-02	14
EuRho	Clobetasol propionate	1408018-01	16
EuRho	Clotrimazole	1108053-02	2
EuRho	Dang Shen	g077h2001222	2
EuRho	Dithranol	1203025-01	1
EuRho	Fructose	171003401	3
EuRho	Gentamicin sulfate	1808013-01	5
EuRho	Gentamicin sulfate	FG1112482	1
EuRho	Glucose monohydrate	1108061-02	1
EuRho	Glucose monohydrate	2710Q-02209	1
EuRho	Glucose monohydrate	151780	2
EuRho	Hydrocortisone	1405001-01	6
EuRho	Hydroxyethylcellulose...	3081202	1
EuRho	Lactose monohydrate	1112020-01	6
EuRho	Lactose monohydrate	170301901	3
EuRho	Lidocaine hydrochlori...	161201502	4
EuRho	Lidocaine hydrochlori...	1108017-01	6
EuRho	Lidocaine hydrochlori...	140302501	3
EuRho	Macrogols 4000	1206009-01	1
EuRho	Magnesium sulfate hep...	3512M-01835	1
EuRho	Magnesium sulfate hep...	1303008-01	5
EuRho	Mannitol	1502022-02	6
EuRho	Mannitol	1911M-01844	1
EuRho	Menthol	2910A-01855	1
EuRho	Menthol	915A-01855	1
EuRho	Methyl-4-hydroxybenzo...	110403101	1
EuRho	Metronidazole	1406010-01	21
EuRho	Metronidazole	1212001-02	18
EuRho	Metronidazole	110903002	1
EuRho	Minoxidil	1105030-02	1
EuRho	Minoxidil	1108039-02	2
EuRho	Naphazoline hydrochlo...	1114I-001885	1
EuRho	Prasterone	1311002-01	1
EuRho	Prednisolone acetate	1306040-03	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
EuRho	Prednisolone, microni...	1012012-02	4
EuRho	Prednisolone, microni...	1012012-01	1
EuRho	Riboflavin	1810005-01	1
EuRho	Salicylic acid	Retest:13052403	1
EuRho	Salicylic acid	1206007-01	10
EuRho	Sodium citrate	1070102901	1
EuRho	Sodium citrate	3032981	2
EuRho	Sodium citrate	1107011-02	3
EuRho	Sodium citrate	1211030-01	1
EuRho	Sodium citrate	1205026-02	5
EuRho	Sodium citrate	4611M-01900	1
EuRho	Sodium citrate	1412014-01	2
EuRho	Sodium sulfate decahy...	1210005-01	2
EuRho	Sorbitol	161001802	1
EuRho	Talcum	1105024-02	1
EuRho	Tetracaine hydrochlor...	1615I-02039	1
EuRho	Tetracaine hydrochlor...	2906009V-04	1
EuRho	Tetracaine hydrochlor...	1110010-01	8
EuRho	Tetracycline hydrochl...	1402030-02	10
EuRho	Tosylchloramide sodium	133610607	1
EuRho	Tretinoin	151007-01	1
EuRho	Triamcinolone acetonide	1508008-02	18
EuRho	Triamcinolone acetonide	150800-08-02	1
EuRho	Triamcinolone acetonide	NW2518M	2
EuRho	Triclosan	151004501	3
EuRho	Triclosan	1107030-01	5
EuRho	Triclosan	1202018-01	2
EuRho / Noweda	Minoxidil	1304029-02	20
EuRho / Phoenix	Allantoin	1406011-03	1
EuRho / Phoenix	Sodium citrate	151101017-01	1
EuRho Anzag ek:3,47EUR	Cetyl alcohol	110505501	1
EuRho Pharm	Aluminium chloride he...	1109025-02	3
EuRho/ Anzag	Dequalinium chloride	80(1804009-01)	1
EuRho/AHD	Tetracycline hydrochl...	1110029-01/Nachprüfu	1
EuRho/Anzag	Sorbitol	4712I-02298	1
EuRho/Jenne	Diclofenac sodium	1(1711026-03)	1
EuRho/Kehr	Minoxidil	30/17	1
EuRho/Noweda	Metronidazole	180601M1706022-01	1
Eurho/Phoenix	Naphazoline hydrochlo...	4413I-01885	1
EuRho/Phoenix	Hydroxyethylcellulose...	1804003-01	7
EuRho/Phönix	Chloramphenicol	123012-01	1
EuRho/Phönix	Clotrimazole	1210039-2	1
EuRho/Phönix	Erythromycin	1205035	1
EuRho/Phönix	Magnesium hydrogen ci...	1202034-02	1
EuRho/Phönix	Mometasone furoate, m...	1307037-01	3
Eurho/Sanacorp	Glucose monohydrate	1108061-01	6
EuRho/Spangro	Diclofenac sodium	1711026-02	7
Euro	Clioquinol	1509038-01	13
Euro	Salicylic acid	1509007-02	15
Euro	Sodium citrate	1815M-01900	1
Euro	Tretinoin	14111016-01	1
euro AHD	Betamethasone dipropi...	1801045-015787	1
euro AHD	Betamethasone dipropi...	1801045-015715	1
euro AHD	Clobetasol propionate	1804002-015788	1
euro AHD	Clobetasol propionate	1704020-015592	1
euro AHD	Gentamicin sulfate	1711025-015586	1
euro AHD	Miconazole nitrate	1710020-015639	1
euro AHD	Miconazole nitrate	1710020-015633	1
euro AHD	Salicylic acid	1707009-025709	1
Euro / Jenne	Magnesium sulfate hep...	1401029-01	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
euro / Phoenix 19.02....	Citric acid	11609021	1
EuRo Gehe	Acetylsalicylic acid	147026-01	1
Euro OCT	Sodium hydrogen carbo...	1307017-01	1
Euro OCT/Noweda	Tobramycin sulfate	1510035-01	1
Euro OT	Clobetasol propionate	14081018-02	1
Euro OT/Phönix	Potassium citrate	3036325	1
Euro OT/Phönix	Tetracaine hydrochlor...	TH0030315	2
Euro OTC	Acetylcysteine	120802601	2
Euro OTC	Acetylcysteine	1510016-01	2
Euro OTC	Acetylcysteine	1505011-01	1
Euro OTC	Acetylcysteine	12299010	1
Euro OTC	Acetylsalicylic acid	1601042-02	1
Euro OTC	Acetylsalicylic acid	1011015-01	1
Euro OTC	Acetylsalicylic acid	1107007-01	1
Euro OTC	Acetylsalicylic acid	1212005-01	4
Euro OTC	Acetylsalicylic acid	1407025-01	17
Euro OTC	Acetylsalicylic acid	1407025	1
Euro OTC	Acetylsalicylic acid	1307024-01	9
Euro OTC	Acetylsalicylic acid	1501046-01	4
Euro OTC	Acetylsalicylic acid	20030611-1	1
Euro OTC	Acetylsalicylic acid	815A-01012	1
Euro OTC	Alfatradiol	1404001-02	1
Euro OTC	Alfatradiol	AP33-16-1-62P2985A	1
Euro OTC	Alfatradiol	160301901	1
Euro OTC	Alfatradiol	1404001-01	2
Euro OTC	Alfatradiol	121696	1
Euro OTC	Allantoin	1401046-01	4
Euro OTC	Allantoin	1406011-01	7
Euro OTC	Allantoin	1406011-05	2
Euro OTC	Allantoin	1406011-02	1
Euro OTC	Allantoin	1609012-03	1
Euro OTC	Allantoin	160901202	1
Euro OTC	Allantoin	1406011-04	2
Euro OTC	Allantoin	1406011	2
Euro OTC	Allantoin	1609012-01	8
Euro OTC	Allantoin	140601103	1
Euro OTC	Allantoin	140104601	1
Euro OTC	Almonds, bitter, whole	290973	2
Euro OTC	alpha-tocopheryl acet...	181002401	1
Euro OTC	alpha-tocopheryl acet...	1705003-04	1
Euro OTC	alpha-tocopheryl acet...	170500301	1
Euro OTC	alpha-tocopheryl acet...	1110021-02	2
Euro OTC	alpha-tocopheryl acet...	1603046-02	1
Euro OTC	alpha-tocopheryl acet...	1507001-01	2
Euro OTC	alpha-tocopheryl acet...	160304602	1
Euro OTC	Aluminium chloride he...	1103007-01	1
Euro OTC	Aluminium chloride he...	1109025-01	1
Euro OTC	Aluminium chloride he...	1507002-01	1
Euro OTC	Aluminium chloride he...	co31a	1
Euro OTC	Aluminium chloride he...	1406021-01	1
Euro OTC	Aluminium chloride he...	120402602	1
Euro OTC	Aluminium chloride he...	1210014-02	29
Euro OTC	Aluminium chloride he...	1204026-01	20
Euro OTC	Aluminium chloride he...	1210014-01	25
Euro OTC	Aluminium chloride he...	1213M-01095	1
Euro OTC	Aluminium chloride he...	1307011-01	10
Euro OTC	Aluminium chloride he...	120402601	3
Euro OTC	Aluminium chloride he...	121001402	1
Euro OTC	Aluminium chloride he...	2111Q-01095	1
Euro OTC	Aluminium chloride he...	1307011	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Aluminium chloride he...	1210012-01	1
Euro OTC	Aluminium chloride he...	1307011-02	2
Euro OTC	Aluminium chloride he...	1210014-02-336331909	1
Euro OTC	Aluminium potassium s...	1004046-02	1
Euro OTC	Ambroxol hydrochloride	1711021-01	12
Euro OTC	Ambroxol hydrochloride	1702004-01	7
Euro OTC	Ambroxol hydrochloride	1702004-02	9
Euro OTC	Ambroxol hydrochloride	27061702	1
Euro OTC	Ambroxol hydrochloride	160215-02	1
Euro OTC	Ambroxol hydrochloride	170200402	1
Euro OTC	Ambroxol hydrochloride	171102101	1
Euro OTC	Ambroxol hydrochloride	1806009	1
Euro OTC	Ambroxol hydrochloride	1711021-02	17
Euro OTC	Ambroxol hydrochloride	1806009-01	13
Euro OTC	Ambroxol hydrochloride	AXL0592201803	1
Euro OTC	Ambroxol hydrochloride	1(1806009-01)	1
Euro OTC	Ambroxol hydrochloride	181103901	3
Euro OTC	Ambroxol hydrochloride	1811039-01	7
Euro OTC	Ambroxol hydrochloride	1904001-01	1
Euro OTC	Ambroxol hydrochloride	171102-01	1
Euro OTC	Ambroxol hydrochloride	1604032-01	1
Euro OTC	Ambroxol hydrochloride	1609002-01	8
Euro OTC	Ambroxol hydrochloride	ax10492201608	1
Euro OTC	Ambroxol hydrochloride	1604035-01	1
Euro OTC	Ambroxol hydrochloride	1311032-01	1
Euro OTC	Ambroxol hydrochloride	2(1904001-01)	1
Euro OTC	Aminomix NAC (H1002)	16001017/0	2
Euro OTC	Amitriptyline hydroch...		1
Euro OTC	Amitriptyline hydroch...	2715E-07209	1
Euro OTC	Ammonium chloride	1303029-01	1
Euro OTC	Ammonium chloride	1203029-01	1
Euro OTC	Ammonium chloride	1212017-01	2
Euro OTC	Ammonium chloride	512A-01104	1
Euro OTC	Amphotericin B	1305033-03	12
Euro OTC	Angelica root	301767	2
Euro OTC	Antioxidants NAT comp...	16000970	3
Euro OTC	Arginine	L1206006-01	1
Euro OTC	Arginine	1604003-01	1
Euro OTC	Arginine	1902024-01	1
Euro OTC	Arginine	1710010-01	2
Euro OTC	Arginine	1105040-01	1
Euro OTC	Arginine hydrochloride	L1506023-01	1
Euro OTC	Arginine hydrochloride	L1407040-01	1
Euro OTC	Arginine hydrochloride	L1212009-01	3
Euro OTC	Arginine hydrochloride	1306023-01	2
Euro OTC	Arginine hydrochloride	L1407040-02	1
Euro OTC	Ascorbic acid	12061408	1
Euro OTC	Ascorbic acid	17041406	1
Euro OTC	Ascorbic acid	L1202011-02	1
Euro OTC	Ascorbic acid	17011404	1
Euro OTC	Ascorbic acid	1710032-01	2
Euro OTC	Ascorbic acid	10071418	1
Euro OTC	Ascorbic acid	11121302	1
Euro OTC	Ascorbic acid	L60502801	1
Euro OTC	Ascorbic acid	28081303	1
Euro OTC	Ascorbyl palmitate	16I30-B03	1
Euro OTC	Bacitracin	17L11-B02-348461	1
Euro OTC	Beclometasone dipropi...	1502032-02	2
Euro OTC	Beclometasone dipropi...	C0783111020	1
Euro OTC	Beclometasone dipropi...	1203022-03	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Beclometasone dipropi...	1203022-01	2
Euro OTC	Benzocaine	15227313	2
Euro OTC	Benzocaine	13L06-B02-297505	1
Euro OTC	Betaine hydrochloride	1202008-01	1
Euro OTC	Betaine hydrochloride	1305027-01	2
Euro OTC	Betaine hydrochloride	L1412019-01	1
Euro OTC	Betaine hydrochloride	1510011-02	1
Euro OTC	Betaine hydrochloride	L1701040-01	1
Euro OTC	Betaine hydrochloride	L1510011-03	2
Euro OTC	Betamethasone dipropi...	1801045-01	9
Euro OTC	Betamethasone dipropi...	160103401	1
Euro OTC	Betamethasone dipropi...	1610014-015308	1
Euro OTC	Betamethasone dipropi...	1311023-01	2
Euro OTC	Betamethasone dipropi...	13372501	4
Euro OTC	Betamethasone dipropi...	131102301	1
Euro OTC	Betamethasone dipropi...	180808185	1
Euro OTC	Betamethasone dipropi...	1610014-01	24
Euro OTC	Betamethasone dipropi...	150300302	1
Euro OTC	Betamethasone dipropi...	121002001	1
Euro OTC	Betamethasone dipropi...	1303006-01	1
Euro OTC	Betamethasone dipropi...	1503003-02	9
Euro OTC	Betamethasone dipropi...	1503003-01	10
Euro OTC	Betamethasone dipropi...	1311023-02	4
Euro OTC	Betamethasone dipropi...	1506016-02	1
Euro OTC	Betamethasone dipropi...	1311023-01	6
Euro OTC	Betamethasone dipropi...	17B16-B04-332218	1
Euro OTC	Betamethasone dipropi...	131102301	1
Euro OTC	Betamethasone dipropi...	1503003-01	1
Euro OTC	Betamethasone dipropi...	1305041-01	1
Euro OTC	Betamethasone dipropi...	1601034-01	21
Euro OTC	Betamethasone dipropi...	13125702	1
Euro OTC	Betamethasone dipropi...	1303006/01	1
Euro OTC	Betamethasone dipropi...	130541-01	1
Euro OTC	Betamethasone dipropi...	1305041-01	8
Euro OTC	Betamethasone dipropi...	180104501	2
Euro OTC	Betamethasone dipropi...	1303006-02	1
Euro OTC	Betamethasone dipropi...	170502001	2
Euro OTC	Betamethasone dipropi...	1303006-01	6
Euro OTC	Betamethasone dipropi...	1305041-02	3
Euro OTC	Betamethasone valerate	2014Q1181	1
Euro OTC	Betamethasone valerate	1311039-01	22
Euro OTC	Betamethasone valerate	1704018-01	1
Euro OTC	Betamethasone valerate	1704019-0	1
Euro OTC	Betamethasone valerate	16270505	15
Euro OTC	Betamethasone valerate	17K16-B08-345541	1
Euro OTC	Betamethasone valerate	1108045-01	2
Euro OTC	Betamethasone valerate	161001001	1
Euro OTC	Betamethasone valerate	BV116/1117	1
Euro OTC	Betamethasone valerate	120801404	1
Euro OTC	Betamethasone valerate	131103901	1
Euro OTC	Betamethasone valerate	140601903	3
Euro OTC	Betamethasone valerate	1704019-01	95
Euro OTC	Betamethasone valerate	1311039	1
Euro OTC	Betamethasone valerate	1406019-01	27
Euro OTC	Betamethasone valerate	1704019-01	1
Euro OTC	Betamethasone valerate	1206014-02	1
Euro OTC	Betamethasone valerate	14C19-B02-292650	5
Euro OTC	Betamethasone valerate	131139-01	1
Euro OTC	Betamethasone valerate	1704019	2
Euro OTC	Betamethasone valerate	1704019-1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Betamethasone valerate	1801024-01	62
Euro OTC	Betamethasone valerate	140601901	1
Euro OTC	Betamethasone valerate	1406011903	1
Euro OTC	Betamethasone valerate	160623-01	1
Euro OTC	Betamethasone valerate	16100100-01	1
Euro OTC	Betamethasone valerate	120812-02	1
Euro OTC	Betamethasone valerate	1208014-01	16
Euro OTC	Betamethasone valerate	3313U-01181	1
Euro OTC	Betamethasone valerate	1510037	1
Euro OTC	Betamethasone valerate	1606023-01	39
Euro OTC	Betamethasone valerate	151003701	1
Euro OTC	Betamethasone valerate	110804503	1
Euro OTC	Betamethasone valerate	1406019-03	37
Euro OTC	Betamethasone valerate	120801403	1
Euro OTC	Betamethasone valerate	1208014-03	5
Euro OTC	Betamethasone valerate	BV-BLT221/0314	1
Euro OTC	Betamethasone valerate	1811008-1	1
Euro OTC	Betamethasone valerate	1811008-01	4
Euro OTC	Betamethasone valerate	110804-02	1
Euro OTC	Betamethasone valerate	140601902	1
Euro OTC	Betamethasone valerate	2507211181	1
Euro OTC	Betamethasone valerate	170419-01	1
Euro OTC	Betamethasone valerate	1208014-02	14
Euro OTC	Betamethasone valerate	1108045-02	3
Euro OTC	Betamethasone valerate	13K29-B32-289364	1
Euro OTC	Betamethasone valerate	1406019-02	26
Euro OTC	Bifonazole	1206027-01	2
Euro OTC	Bifonazole	1206027-02	1
Euro OTC	Bifonazole	1510001-02	1
Euro OTC	Bifonazole	170304101	1
Euro OTC	Bifonazole	13G10N02	1
Euro OTC	Bifonazole	16703-B08-333937	1
Euro OTC	Bifonazole	1510001-02	4
Euro OTC	Biotin	24011801	1
Euro OTC	Biotin	5121812	1
Euro OTC	Biotin	1111039-02	1
Euro OTC	Biotin	1301026-01	3
Euro OTC	Biotin	1802022-01	3
Euro OTC	Biotin	3021501	1
Euro OTC	Biotin	2041517	1
Euro OTC	Biotin	28071506	1
Euro OTC	Biotin	1411047-01	1
Euro OTC	Biotin	1603017-01	1
Euro OTC	Biotin	161000301	1
Euro OTC	Biotin	1610003-01	4
Euro OTC	Biotin	10061709	1
Euro OTC	Biotin	1402028-01	4
Euro OTC	Biotin	1708007-01	1
Euro OTC	Budesonide, micronized	1604018-01	1
Euro OTC	Budesonide, micronized	1311022-01	4
Euro OTC	Burdock root, cut	13122803	1
Euro OTC	Caffeine	291001903	2
Euro OTC	Caffeine	1411015-01	3
Euro OTC	Caffeine	325241AX20	1
Euro OTC	Caffeine	1703006-01	19
Euro OTC	Caffeine	180881701	1
Euro OTC	Caffeine	1603009-02	4
Euro OTC	Caffeine	1211012-01	1
Euro OTC	Caffeine	1211012-02	1
Euro OTC	Caffeine	1603009-01	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Caffeine	160300901	2
Euro OTC	Caffeine	1603009/01	1
Euro OTC	Caffeine	1808017-01	3
Euro OTC	Caffeine	1411015-01	1
Euro OTC	Caffeine	1503020-02	8
Euro OTC	Caffeine	1703006-02	2
Euro OTC	Caffeine	1307049-01	3
Euro OTC	Caffeine	26081406	1
Euro OTC	Caffeine	16462119	1
Euro OTC	Caffeine	1109017-01	1
Euro OTC	Calcium ascorbate	1709004-01	1
Euro OTC	Calcium ascorbate	1709005	1
Euro OTC	Calcium ascorbate	1801010	1
Euro OTC	Calcium citrate tetra...	3037086	1
Euro OTC	Calcium citrate tetra...	1606001-02	2
Euro OTC	Calcium citrate tetra...	140801701	1
Euro OTC	Calcium citrate tetra...	1606001-01	1
Euro OTC	Calcium citrate tetra...	1502021-02	1
Euro OTC	Calcium citrate tetra...	1211007-01	4
Euro OTC	Calcium citrate tetra...	3006222	1
Euro OTC	Calcium citrate tetra...	1707025-01	1
Euro OTC	Calcium citrate tetra...	1502012-01	1
Euro OTC	Calcium citrate tetra...	1009030-01	1
Euro OTC	Calcium citrate tetra...	1701001-01	1
Euro OTC	Calcium citrate tetra...	1509039-02	7
Euro OTC	Calcium citrate tetra...	1803031-02	1
Euro OTC	Calcium citrate tetra...	1058904	1
Euro OTC	Calcium citrate tetra...	1207006-01	3
Euro OTC	Calcium citrate tetra...	15177704	1
Euro OTC	Calcium citrate tetra...	1402003-01	3
Euro OTC	Calcium citrate tetra...	1408017-01	5
Euro OTC	Calcium citrate tetra...	13031204	1
Euro OTC	Calcium citrate tetra...	1305018-02	1
Euro OTC	Calcium citrate tetra...	1701001-02	1
Euro OTC	Calcium citrate tetra...	1305018-01	4
Euro OTC	Calcium citrate tetra...	1710007-01	1
Euro OTC	Calcium gluconate	115027-01	1
Euro OTC	Calcium gluconate	140704102	2
Euro OTC	Calcium gluconate	140704101	1
Euro OTC	Calcium gluconate	1606034-02	6
Euro OTC	Calcium gluconate	1407041-01	1
Euro OTC	Calcium gluconate	1304012-02	2
Euro OTC	Calcium gluconate	1210001287	1
Euro OTC	Calcium gluconate	1407041-02	10
Euro OTC	Calcium gluconate	1304012-01	6
Euro OTC	Calcium gluconate	1209006-01	2
Euro OTC	Calcium gluconate	1304012-01	2
Euro OTC	Calcium gluconate	20270912-1	1
Euro OTC	Calcium gluconate	1602001644	1
Euro OTC	Calcium gluconate	513-Q1211	1
Euro OTC	Calcium gluconate	1712007-01	4
Euro OTC	Calcium gluconate	1105027-01	1
Euro OTC	Calcium gluconate	1105027-02	3
Euro OTC	Calcium gluconate	130401202	1
Euro OTC	Calcium gluconate	1209006-01	9
Euro OTC	Calcium orotate dihyd...	1401006-01	1
Euro OTC	Calcium orotate dihyd...	1203037-01	2
Euro OTC	Calcium pantothenate	1711003-01	4
Euro OTC	Calcium pantothenate	1209022-02pzn	1
Euro OTC	Calcium pantothenate	1111025-01	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Calcium pantothenate	1209022-01	1
Euro OTC	Calcium pantothenate	1305034-01	4
Euro OTC	Calcium pantothenate	1803014-01	3
Euro OTC	Calcium pantothenate	1607023-01	4
Euro OTC	Calcium pantothenate	18111707	1
Euro OTC	Calcium pantothenate	1410011-01	2
Euro OTC	Calcium pantothenate	10021505	1
Euro OTC	Calcium pantothenate	5121405	1
Euro OTC	Calcium pantothenate	1209022-02	1
Euro OTC	Calcium pantothenate	1702018-01	1
Euro OTC	Calcium pantothenate	1312027-01	3
Euro OTC	Calcium pantothenate	1410011-02	4
Euro OTC	Calcium pantothenate	1603042-01	3
Euro OTC	Calcium pantothenate	130503401	2
Euro OTC	Calcium pantothenate	17021401	1
Euro OTC	Camphor	1101029-02	1
Euro OTC	Camphor	1309019-02	1
Euro OTC	Camphor	1202021	1
Euro OTC	Camphor	821509	1
Euro OTC	Camphor	1309019-03	1
Euro OTC	Camphor	1202021-02	2
Euro OTC	Camphor	1406029-01	2
Euro OTC	Camphor	1303013-01	5
Euro OTC	Camphor	1303013-02	11
Euro OTC	Camphor	1210030-01	1
Euro OTC	Camphor	1101029-02/Ch241011	1
Euro OTC	Camphor	113Q-01221	1
Euro OTC	Camphor	1406029-02	2
Euro OTC	Camphor	130313-02	1
Euro OTC	Camphor	2803020-01	1
Euro OTC	Camphor	1210030-02	4
Euro OTC	Camphor	121003001	1
Euro OTC	Camphor	1210030-02	1
Euro OTC	Camphor	1210030-01	7
Euro OTC	Camphor	121003002	1
Euro OTC	Capsaicin	1806013-01	1
Euro OTC	Capsaicin	1703023-025346	1
Euro OTC	Capsaicin	1401043-02	2
Euro OTC	Capsaicin	1603005-01	13
Euro OTC	Captopril	17210204	3
Euro OTC	Captopril	15189507	3
Euro OTC	Carbomers 35000 / 50000	1210037-01	6
Euro OTC	Carbomers 35000 / 50000	1202020-01	2
Euro OTC	Carbomers 35000 / 50000	4311027-01	1
Euro OTC	Carbomers 35000 / 50000	101217279	1
Euro OTC	Carbomers 35000 / 50000	1311027-01	4
Euro OTC	Carbomers 35000 / 50000	1512010-01	8
Euro OTC	Carbomers 35000 / 50000	1410001-01	16
Euro OTC	Carbomers 35000 / 50000	1604032-01	9
Euro OTC	Carbomers 35000 / 50000	12052014A	1
Euro OTC	Carbomers 35000 / 50000	120202001	2
Euro OTC	Carbomers 35000 / 50000	1610024-01	15
Euro OTC	Carbomers 35000 / 50000	1512010-02	4
Euro OTC	Carbomers 35000 / 50000	101616773	2
Euro OTC	Carbomers 35000 / 50000	1706015-01	2
Euro OTC	Carbomers 35000 / 50000	141000101	1
Euro OTC	Carbomers 35000 / 50000	1210037-01	1
Euro OTC	Carbomers 35000 / 50000	141000-01	1
Euro OTC	Carbomers 35000 / 50000	1607028-01	8
Euro OTC	Carbomers 35000 / 50000	700143(Artikelnr.)	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Carbomers 35000 / 50000	151210-02	1
Euro OTC	Carbomers 35000 / 50000	1801044-01	4
Euro OTC	Carbomers 35000 / 50000	102026590	1
Euro OTC	Carbomers 35000 / 50000	41(1810019-01)	1
Euro OTC	Carbomers 35000 / 50000	130800601	1
Euro OTC	Carbomers 35000 / 50000	1810019-01	1
Euro OTC	Carbomers 35000 / 50000	1311027-02	8
Euro OTC	Carbomers 35000 / 50000	1305007-01	1
Euro OTC	Carbomers 35000 / 50000	1702012-01	1
Euro OTC	Cellulose, microcryst...	1712009-01	1
Euro OTC	Cellulose, microcryst...	1610030-01	3
Euro OTC	Cellulose, microcryst...	150722-01	1
Euro OTC	Cellulose, microcryst...	135010001	1
Euro OTC	Cellulose, microcryst...	1607042-01	2
Euro OTC	Cellulose, microcryst...	1604029-01	7
Euro OTC	Cellulose, microcryst...	150702201	1
Euro OTC	Cellulose, microcryst...	1410002-01	4
Euro OTC	Cellulose, microcryst...	121103701	1
Euro OTC	Cellulose, microcryst...	700204	1
Euro OTC	Cellulose, microcryst...	1309001-01	3
Euro OTC	Cellulose, microcryst...	1211037-01	3
Euro OTC	Cellulose, microcryst...	160704201	1
Euro OTC	Cellulose, microcryst...	1712009-02	2
Euro OTC	Cetostearyl alcohol	1107035-01	1
Euro OTC	Cetostearyl alcohol	1504010-01	1
Euro OTC	Cetostearyl alcohol	1404019-01	1
Euro OTC	Cetostearyl alcohol	1204015-01	1
Euro OTC	Cetostearyl alcohol	15271403	1
Euro OTC	Cetostearyl alcohol	1305035-01	1
Euro OTC	Cetostearyl alcohol	1404019-02	1
Euro OTC	Cetostearyl alcohol	20270912-2	1
Euro OTC	Cetostearyl alcohol	17050001-01	1
Euro OTC	Cetostearyl alcohol	1604013-02	3
Euro OTC	Cetostearyl alcohol	1510043-01	1
Euro OTC	Cetostearyl alcohol	150401001	1
Euro OTC	Cetostearyl alcohol	1705001-01	2
Euro OTC	Cetyl alcohol	1206020-01	2
Euro OTC	Cetyl alcohol	1211002-01	2
Euro OTC	Cetyl alcohol	150818C	1
Euro OTC	Cetyl alcohol	1105055-01	1
Euro OTC	Cetyl alcohol	1508014-01	1
Euro OTC	Cetyl palmitate 15	160419C	1
Euro OTC	Chloral hydrate	16122001	1
Euro OTC	Chloramphenicol	1704021-01,20180830M	1
Euro OTC	Chloramphenicol	1407043-01	22
Euro OTC	Chloramphenicol	1704021-01	67
Euro OTC	Chloramphenicol	314Q1259	1
Euro OTC	Chloramphenicol	1704021-02	28
Euro OTC	Chloramphenicol	1903028-01	8
Euro OTC	Chloramphenicol	1602002-04	1
Euro OTC	Chloramphenicol	170402103	1
Euro OTC	Chloramphenicol	1704021	1
Euro OTC	Chloramphenicol	1308017.01	1
Euro OTC	Chloramphenicol	1704021-04	14
Euro OTC	Chloramphenicol	170402101	1
Euro OTC	Chloramphenicol	130801701	1
Euro OTC	Chloramphenicol	1704021-015505	1
Euro OTC	Chloramphenicol	T02-55/1407043-01	1
Euro OTC	Chloramphenicol	1903028-02	2
Euro OTC	Chloramphenicol	160200203	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Chloramphenicol	150203602	1
Euro OTC	Chloramphenicol	1308017-01	13
Euro OTC	Chloramphenicol	1509211259	1
Euro OTC	Chloramphenicol	B:1006047-02	1
Euro OTC	Chloramphenicol	1203012-01	7
Euro OTC	Chloramphenicol	DY00914151	2
Euro OTC	Chloramphenicol	110700202	1
Euro OTC	Chloramphenicol	1502036022016030911	1
Euro OTC	Chloramphenicol	1107002-01	2
Euro OTC	Chloramphenicol	1602002-03	27
Euro OTC	Chloramphenicol	1012013-01	2
Euro OTC	Chloramphenicol	1407043	1
Euro OTC	Chloramphenicol	1602002-02	27
Euro OTC	Chloramphenicol	1502036	1
Euro OTC	Chloramphenicol	1502036-01	38
Euro OTC	Chloramphenicol	1602002-01	40
Euro OTC	Chloramphenicol	150203601	5
Euro OTC	Chloramphenicol	150203603	2
Euro OTC	Chloramphenicol	1502036-03	21
Euro OTC	Chloramphenicol	1502036-02	15
Euro OTC	Chloramphenicol	1602002-07	1
Euro OTC	Chloramphenicol	1203012-12	1
Euro OTC	Chloramphenicol	130801702	1
Euro OTC	Chloramphenicol	16020002-03	1
Euro OTC	Chloramphenicol	1209019-02	14
Euro OTC	Chloramphenicol	1602002-03,20170427M	1
Euro OTC	Chloramphenicol	120901902	3
Euro OTC	Chloramphenicol	1209010-02	3
Euro OTC	Chloramphenicol	1302024-01	14
Euro OTC	Chloramphenicol	1209019-01	6
Euro OTC	Chloramphenicol	1203012-02	12
Euro OTC	Chloramphenicol	1302024-02	8
Euro OTC	Chloramphenicol	2121211203012-0	1
Euro OTC	Chlorhexidine diacetate	10154430	1
Euro OTC	Chlorhexidine diacetate	15261004	1
Euro OTC	Chlorhexidine diacetate	1785805	1
Euro OTC	Chlortetracycline hyd...	10I29-N04	1
Euro OTC	Chlortetracycline hyd...	16b10-b09-326865	1
Euro OTC	Chlortetracycline hyd...	1412008/03	1
Euro OTC	Chlortetracycline hyd...	2901009V-03	1
Euro OTC	Chlortetracycline hyd...	1803037-01	1
Euro OTC	Chlortetracycline hyd...	1108051-01	9
Euro OTC	Chlortetracycline hyd...	110805102	2
Euro OTC	Chlortetracycline hyd...	5113I01155	1
Euro OTC	Chlortetracycline hyd...	4612Q-01155	1
Euro OTC	Chlortetracycline hyd...	1412008-01	10
Euro OTC	Chlortetracycline hyd...	1412008-01	1
Euro OTC	Chlortetracycline hyd...	1412008-03	1
Euro OTC	Chlortetracycline hyd...	1412008-02	16
Euro OTC	Chlortetracycline hyd...	1412008-03	4
Euro OTC	Chlortetracycline hyd...	141200802	3
Euro OTC	Chlortetracycline hyd...	1412008-02	1
Euro OTC	Chlortetracycline hyd...	1604015-01	30
Euro OTC	Chlortetracycline hyd...	160401502	1
Euro OTC	Cholecalciferol dry c...	1806003-01	28
Euro OTC	Cholecalciferol dry c...	1707027-01	17
Euro OTC	Cholecalciferol dry c...	1508032-01	2
Euro OTC	Cholecalciferol dry c...	1012025-01	1
Euro OTC	Cholecalciferol dry c...	1602019-01	2
Euro OTC	Cholecalciferol dry c...	1502021-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Cholecalciferol dry c...	170153-01	5
Euro OTC	Cholecalciferol dry c...	L1403012-01	1
Euro OTC	Cholecalciferol dry c...	1802014-01	16
Euro OTC	Cholecalciferol dry c...	150803201	1
Euro OTC	Cholecalciferol dry c...	1210051-01	1
Euro OTC	Cholecalciferol dry c...	1701053-01	11
Euro OTC	Cholecalciferol dry c...	1801038-01	4
Euro OTC	Choline hydrogen tart...	1309023-01	2
Euro OTC	Choline hydrogen tart...	1711022-01	1
Euro OTC	Choline hydrogen tart...	1707023-01	1
Euro OTC	Choline hydrogen tart...	1607020-01	2
Euro OTC	Choline hydrogen tart...	1603044-01	1
Euro OTC	Chromium yeast 0.2%	L1311026-01	1
Euro OTC	Ciclopirox olamine	1718806	1
Euro OTC	Citric acid	L1604061-02	3
Euro OTC	Citric acid	L130404201	2
Euro OTC	Citric acid	L1408011-02	2
Euro OTC	Citric acid	L1508013-02	1
Euro OTC	Citric acid	L1304042.02	1
Euro OTC	Citric acid	L1609021-02	4
Euro OTC	Citric acid	L1304042-02	1
Euro OTC	Citric acid	15.01.15	2
Euro OTC	Citric acid	1505029B	1
Euro OTC	Citric acid	1406014-02	1
Euro OTC	Citric acid	L1408011-01	1
Euro OTC	Citric acid	13396515	3
Euro OTC	Citric acid	702260	3
Euro OTC	Citric acid	L1508013-01	5
Euro OTC	Citric acid	L1609021-01	7
Euro OTC	Citric acid	13175410	5
Euro OTC	Citric acid	L1604062-01	1
Euro OTC	Citric acid	1011047-03	2
Euro OTC	Citrulline	16040038-01	1
Euro OTC	Clindamycin hydrochlo...	1111058-01	7
Euro OTC	Clindamycin hydrochlo...	2908019V-03	5
Euro OTC	Clindamycin hydrochlo...	171001301	1
Euro OTC	Clindamycin hydrochlo...	4988149	1
Euro OTC	Clindamycin hydrochlo...	1401401-01	1
Euro OTC	Clindamycin hydrochlo...	111158-01	1
Euro OTC	Clindamycin hydrochlo...	1111058-02	1
Euro OTC	Clindamycin hydrochlo...	T03-05/1401041-01	1
Euro OTC	Clindamycin hydrochlo...	1111058-03	4
Euro OTC	Clindamycin hydrochlo...	1510038-01	18
Euro OTC	Clindamycin hydrochlo...	1401041-02	6
Euro OTC	Clindamycin hydrochlo...	151100401	2
Euro OTC	Clindamycin hydrochlo...	1510038-02	3
Euro OTC	Clindamycin hydrochlo...	700168	1
Euro OTC	Clindamycin hydrochlo...	1511004-01	15
Euro OTC	Clindamycin hydrochlo...	140104102	4
Euro OTC	Clindamycin hydrochlo...	SCHE1107035	1
Euro OTC	Clindamycin hydrochlo...	130709X8	2
Euro OTC	Clindamycin hydrochlo...	16D05-B08-328259	2
Euro OTC	Clindamycin hydrochlo...	140104101	1
Euro OTC	Clindamycin hydrochlo...	1710013-01	17
Euro OTC	Clindamycin hydrochlo...	1511004-002	1
Euro OTC	Clindamycin hydrochlo...	813M-03534	1
Euro OTC	Clindamycin hydrochlo...	25071401	1
Euro OTC	Clindamycin hydrochlo...	150901B	1
Euro OTC	Clindamycin hydrochlo...	1401041-01	14
Euro OTC	Clindamycin hydrochlo...	8061801	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Clindamycin hydrochlo...	1510030-01	1
Euro OTC	Clioquinol	1509038-31	1
Euro OTC	Clioquinol	1102022-01	2
Euro OTC	Clioquinol	1213A-02562	1
Euro OTC	Clioquinol	1309009-01	9
Euro OTC	Clioquinol	PC86	1
Euro OTC	Clioquinol	1301014	1
Euro OTC	Clioquinol	130900902	1
Euro OTC	Clioquinol	1303009-02	1
Euro OTC	Clioquinol	1301014-01	2
Euro OTC	Clioquinol	13010104-03	1
Euro OTC	Clioquinol	140903101	1
Euro OTC	Clioquinol	150903802	3
Euro OTC	Clioquinol	1409031-01	5
Euro OTC	Clioquinol	1509038-02	29
Euro OTC	Clioquinol	1409031-02	12
Euro OTC	Clioquinol	1509038	2
Euro OTC	Clioquinol	1305070109	1
Euro OTC	Clioquinol	1509038-03	23
Euro OTC	Clioquinol	1301014-02	6
Euro OTC	Clioquinol	1201034-02	13
Euro OTC	Clioquinol	1102022-02	6
Euro OTC	Clioquinol	2507212562	1
Euro OTC	Clioquinol	1509038-04	1
Euro OTC	Clioquinol	161202001	2
Euro OTC	Clioquinol	14090131-02	1
Euro OTC	Clioquinol	2195137	1
Euro OTC	Clioquinol	B:2905031-03	1
Euro OTC	Clioquinol	1612020-01	25
Euro OTC	Clioquinol	1612020_01	1
Euro OTC	Clioquinol	1509038-0	1
Euro OTC	Clioquinol	150903-01	1
Euro OTC	Clioquinol	120103402	1
Euro OTC	Clobetasol propionate	1072015D	1
Euro OTC	Clobetasol propionate	EuroOTC	4
Euro OTC	Clobetasol propionate	1804002-01	42
Euro OTC	Clobetasol propionate	1804002201	1
Euro OTC	Clobetasol propionate	1704020-015498	1
Euro OTC	Clobetasol propionate	13042015A	1
Euro OTC	Clobetasol propionate	1108044-01	4
Euro OTC	Clobetasol propionate	1108044-02	7
Euro OTC	Clobetasol propionate	1408018-02	74
Euro OTC	Clobetasol propionate	1408015-01	1
Euro OTC	Clobetasol propionate	CL0-BLT223/1216	1
Euro OTC	Clobetasol propionate	1211021/02	1
Euro OTC	Clobetasol propionate	1108044-037-16	1
Euro OTC	Clobetasol propionate	1408018	1
Euro OTC	Clobetasol propionate	1105044-04	1
Euro OTC	Clobetasol propionate	1606022-01	31
Euro OTC	Clobetasol propionate	160602201	7
Euro OTC	Clobetasol propionate	80616-1	1
Euro OTC	Clobetasol propionate	121102101	2
Euro OTC	Clobetasol propionate	121K021-01	1
Euro OTC	Clobetasol propionate	1108044-03	17
Euro OTC	Clobetasol propionate	16042014C	1
Euro OTC	Clobetasol propionate	1211021-02	35
Euro OTC	Clobetasol propionate	1211021-03	30
Euro OTC	Clobetasol propionate	CL0119/1217	1
Euro OTC	Clobetasol propionate	140808-02	1
Euro OTC	Clobetasol propionate	12082013B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Clobetasol propionate	4091512	1
Euro OTC	Clobetasol propionate	12110-01	1
Euro OTC	Clobetasol propionate	CL0-BL2322/1012	2
Euro OTC	Clobetasol propionate	2111I-02835	1
Euro OTC	Clobetasol propionate	1601046-01	1
Euro OTC	Clobetasol propionate	140801803	1
Euro OTC	Clobetasol propionate	30092014A	1
Euro OTC	Clobetasol propionate	c1o-bl1t223/126	1
Euro OTC	Clobetasol propionate	171107CP/1704020-01	1
Euro OTC	Clobetasol propionate	1211021-04	20
Euro OTC	Clobetasol propionate	130604002	1
Euro OTC	Clobetasol propionate	110804404	2
Euro OTC	Clobetasol propionate	170402001	8
Euro OTC	Clobetasol propionate	1704020-01	88
Euro OTC	Clobetasol propionate	4113Q-02835	1
Euro OTC	Clobetasol propionate	160602202	5
Euro OTC	Clobetasol propionate	1108044	1
Euro OTC	Clobetasol propionate	1211021-01	20
Euro OTC	Clobetasol propionate	1903014-01	3
Euro OTC	Clobetasol propionate	1606022-02	45
Euro OTC	Clobetasol propionate	1108044-04	17
Euro OTC	Clobetasol propionate	140801801	1
Euro OTC	Clotrimazole	1409027-02Mg10g	1
Euro OTC	Clotrimazole	140902702	1
Euro OTC	Clotrimazole	1409027-03	26
Euro OTC	Clotrimazole	21052015B	1
Euro OTC	Clotrimazole	15042014C	1
Euro OTC	Clotrimazole	1409027	7
Euro OTC	Clotrimazole	120304002	1
Euro OTC	Clotrimazole	1306036-01	47
Euro OTC	Clotrimazole	1210039-02	22
Euro OTC	Clotrimazole	121003901	2
Euro OTC	Clotrimazole	12082013C	1
Euro OTC	Clotrimazole	2313M01276	1
Euro OTC	Clotrimazole	21122013A	1
Euro OTC	Clotrimazole	110501402	1
Euro OTC	Clotrimazole	1210039-01	27
Euro OTC	Clotrimazole	700135	2
Euro OTC	Clotrimazole	1105014-01	1
Euro OTC	Clotrimazole	1203040-03	13
Euro OTC	Clotrimazole	1203040-02	11
Euro OTC	Clotrimazole	1209211276	1
Euro OTC	Clotrimazole	3912E-01276	1
Euro OTC	Clotrimazole	1409027-02	43
Euro OTC	Clotrimazole	1306306-01	1
Euro OTC	Clotrimazole	1314Q1276	1
Euro OTC	Clotrimazole	140927-01	1
Euro OTC	Clotrimazole	1309054-02	1
Euro OTC	Clotrimazole	22122014B	1
Euro OTC	Clotrimazole	1409027-01	28
Euro OTC	Clotrimazole	1309034-02Mg:10g	1
Euro OTC	Clotrimazole	1203040-01	2
Euro OTC	Clotrimazole	130636-01	2
Euro OTC	Clotrimazole	1309034	1
Euro OTC	Clotrimazole	130903401	4
Euro OTC	Clotrimazole	15072014A	1
Euro OTC	Clotrimazole	20130512	1
Euro OTC	Clotrimazole	1309034-02	25
Euro OTC	Clotrimazole	1309034-01	24
Euro OTC	Clotrimazole	1203040	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Clotrimazole	2214A-01276	1
Euro OTC	Clotrimazole	130603601	4
Euro OTC	Clotrimazole	1703028-02	34
Euro OTC	Clotrimazole	150003-02	1
Euro OTC	Clotrimazole	180002801	1
Euro OTC	Clotrimazole	1808028_02	1
Euro OTC	Clotrimazole	16234511	16
Euro OTC	Clotrimazole	1607064-01	1
Euro OTC	Clotrimazole	1607046-02	31
Euro OTC	Clotrimazole	1808028-02	17
Euro OTC	Clotrimazole	1607046-0	1
Euro OTC	Clotrimazole	180802802	1
Euro OTC	Clotrimazole	1902004-01	4
Euro OTC	Clotrimazole	H20160502	1
Euro OTC	Clotrimazole	16074046-01	1
Euro OTC	Clotrimazole	1609046-01	1
Euro OTC	Clotrimazole	1409027-04	28
Euro OTC	Clotrimazole	190200401	1
Euro OTC	Clotrimazole	192004-01	1
Euro OTC	Clotrimazole	160704601	7
Euro OTC	Clotrimazole	1607046-01	61
Euro OTC	Clotrimazole	151000	1
Euro OTC	Clotrimazole	1510003-01	55
Euro OTC	Clotrimazole	151000-01	1
Euro OTC	Clotrimazole	15010003-01	1
Euro OTC	Clotrimazole	1409027_03	1
Euro OTC	Clotrimazole	151003-01	1
Euro OTC	Clotrimazole	1510003-02	20
Euro OTC	Clotrimazole	16B11-B04-320088	6
Euro OTC	Clotrimazole	1605003C	1
Euro OTC	Clotrimazole	1510003_01	1
Euro OTC	Clotrimazole	1808028	1
Euro OTC	Clotrimazole	160802801	1
Euro OTC	Clotrimazole	151000301	1
Euro OTC	Clotrimazole	140927-02	1
Euro OTC	Clotrimazole	17053183	1
Euro OTC	Clotrimazole	20160801	5
Euro OTC	Clotrimazole	170302801	8
Euro OTC	Clotrimazole	170302802	1
Euro OTC	Clotrimazole	17A021	1
Euro OTC	Clotrimazole	1703028-03,180626C	1
Euro OTC	Clotrimazole	170302803	2
Euro OTC	Clotrimazole	170615C/1703028-01	1
Euro OTC	Clotrimazole	180802801	2
Euro OTC	Clotrimazole	1808028-01	36
Euro OTC	Clotrimazole	41(1808028-01)	1
Euro OTC	Clotrimazole	1703028-03	27
Euro OTC	Clotrimazole	1703028-02	1
Euro OTC	Clotrimazole	1703028-01	69
Euro OTC	Colistin sulfate	1507044-01	1
Euro OTC	Cyanocobalamin	5061204	1
Euro OTC	Cysteine hydrochlorid...	1304019-01	1
Euro OTC	Cystine	30041809	1
Euro OTC	Cystine	1207013-01	2
Euro OTC	Cystine	2081301	1
Euro OTC	Cystine	1801016-01	2
Euro OTC	Cystine	1204090	1
Euro OTC	Cystine	3021502	1
Euro OTC	Cystine	26081408	1
Euro OTC	Dandelion herb with r...	11030312	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Deoxycholic acid	1305051-01	1
Euro OTC	Dequalinium chloride	1110002-01	3
Euro OTC	Dequalinium chloride	1405012-01	3
Euro OTC	Dequalinium chloride	1210008-03	4
Euro OTC	Dequalinium chloride	1210008-02	1
Euro OTC	Dequalinium chloride	1210008-01	1
Euro OTC	Dequalinium chloride	1103045-01	1
Euro OTC	Dequalinium chloride	15061204	1
Euro OTC	Dequalinium chloride	81(1804009-01)	1
Euro OTC	Dequalinium chloride	1010022-01	2
Euro OTC	Dequalinium chloride	151203201	2
Euro OTC	Dequalinium chloride	1512032-01	16
Euro OTC	Dequalinium chloride	1405012-02	1
Euro OTC	Dequalinium chloride	1210008-02	1
Euro OTC	Dequalinium chloride	1407029-02	5
Euro OTC	Dequalinium chloride	11091701	1
Euro OTC	Dequalinium chloride	1512032-02	3
Euro OTC	Dequalinium chloride	1407029-01	8
Euro OTC	Dequalinium chloride	1512032-01	1
Euro OTC	Dexamethasone	5213A1327	1
Euro OTC	Dexamethasone	130300703	1
Euro OTC	Dexamethasone	1601052-01	24
Euro OTC	Dexamethasone	DZ11005K	1
Euro OTC	Dexamethasone	160105201	2
Euro OTC	Dexamethasone	120900802	2
Euro OTC	Dexamethasone	1601052-02	40
Euro OTC	Dexamethasone	700107	1
Euro OTC	Dexamethasone	160105202	1
Euro OTC	Dexamethasone	161100101	3
Euro OTC	Dexamethasone	1203032-01	4
Euro OTC	Dexamethasone	1303007	1
Euro OTC	Dexamethasone	1512211327	1
Euro OTC	Dexamethasone	130300701	4
Euro OTC	Dexamethasone	1303007-02	10
Euro OTC	Dexamethasone	1111055-02	2
Euro OTC	Dexamethasone	130300702	2
Euro OTC	Dexamethasone	B: 1111055-01	1
Euro OTC	Dexamethasone	111105501	1
Euro OTC	Dexamethasone	DB003/0213	3
Euro OTC	Dexamethasone	L53613	1
Euro OTC	Dexamethasone	1111055-01	5
Euro OTC	Dexamethasone	1107006-01	1
Euro OTC	Dexamethasone	13030007	1
Euro OTC	Dexamethasone	14D11-B01-294031	2
Euro OTC	Dexamethasone	13030007-01	1
Euro OTC	Dexamethasone	1303007-01	33
Euro OTC	Dexamethasone	1303007-03	33
Euro OTC	Dexamethasone	1209008-02	26
Euro OTC	Dexamethasone	1611001-01	45
Euro OTC	Dexamethasone	1708004-01	51
Euro OTC	Dexamethasone	161101-01	1
Euro OTC	Dexamethasone	1611001	1
Euro OTC	Dexamethasone	180208D/1611001-01	1
Euro OTC	Dexamethasone	1806023-02	2
Euro OTC	Dexamethasone	180602301	2
Euro OTC	Dexamethasone	1303007-04	1
Euro OTC	Dexamethasone	1702004-01	1
Euro OTC	Dexamethasone	T45290	1
Euro OTC	Dexamethasone	1708004	1
Euro OTC	Dexamethasone	1806023-01	17

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Dexamethasone	170800401	1
Euro OTC	Dexamethasone	1708004,01	1
Euro OTC	Dexamethasone acetate	1405024-02	6
Euro OTC	Dexamethasone acetate	1110007-01	4
Euro OTC	Dexamethasone acetate	1110007-02	6
Euro OTC	Dexamethasone acetate	1110007-03	1
Euro OTC	Dexamethasone acetate	140502402	1
Euro OTC	Dexamethasone acetate	1405024-03	13
Euro OTC	Dexamethasone acetate	1405024-04	1
Euro OTC	Dexamethasone acetate	1708011-01	7
Euro OTC	Diclofenac sodium	1607016-01	25
Euro OTC	Diclofenac sodium	1503014-02	4
Euro OTC	Diclofenac sodium	1711026-1	1
Euro OTC	Diclofenac sodium	12060250	1
Euro OTC	Diclofenac sodium	1711026-01	6
Euro OTC	Diclofenac sodium	1403013-02	6
Euro OTC	Diclofenac sodium	1607016-02	11
Euro OTC	Diclofenac sodium	140301301	2
Euro OTC	Diclofenac sodium	1403013-01	5
Euro OTC	Diclofenac sodium	1206026	1
Euro OTC	Diclofenac sodium	1206026-02	3
Euro OTC	Diclofenac sodium	12214220	1
Euro OTC	Diclofenac sodium	1206026-01	7
Euro OTC	Diclofenac sodium	15169606	1
Euro OTC	Diclofenac sodium	1607016-03	21
Euro OTC	Diclofenac sodium	13111203	1
Euro OTC	Diclofenac sodium	1111040-01	2
Euro OTC	Diphenylcyclopropenone	17B-13B02-334769	1
Euro OTC	Dithranol	1008011-01	1
Euro OTC	Dithranol	1008011-02	2
Euro OTC	Dithranol	DT/370/2012	1
Euro OTC	Dithranol	B:1008011-01	1
Euro OTC	Dithranol	1203025-02	16
Euro OTC	Dithranol	56/11072013	1
Euro OTC	Emulsifying cetostear...	1510043-01	2
Euro OTC	Emulsifying cetostear...	1305035-01	1
Euro OTC	Emulsifying cetostear...	1404019-03	1
Euro OTC	Emulsifying cetostear...	1604013-02	1
Euro OTC	Emulsifying cetostear...	150410-01	1
Euro OTC	Emulsifying cetostear...	1504010-01	3
Euro OTC	Emulsifying cetostear...	1404019-02	1
Euro OTC	Emulsifying cetostear...	1711001-02	1
Euro OTC	Emulsifying cetostear...	1212019-01	1
Euro OTC	Emulsifying cetostear...	1705001-01	4
Euro OTC	Emulsifying cetostear...	1404019-01	1
Euro OTC	Emulsifying cetostear...	10836573	1
Euro OTC	Equisetum stem, cut	3831A130610	5
Euro OTC	Erythromycin	97(1811019-02)	1
Euro OTC	Erythromycin	1607038-02	67
Euro OTC	Erythromycin	1409019-01	34
Euro OTC	Erythromycin	167038-01	1
Euro OTC	Erythromycin	1409019-02Gehe	1
Euro OTC	Erythromycin	1607038	1
Euro OTC	Erythromycin	1607038-01	63
Euro OTC	Erythromycin	140901901	3
Euro OTC	Erythromycin	15042810	5
Euro OTC	Erythromycin	1508034-02Mg100g	1
Euro OTC	Erythromycin	1409019-01Mg100g	1
Euro OTC	Erythromycin	80616-4	1
Euro OTC	Erythromycin	1508034-02	16

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Erythromycin	1811019-02	2
Euro OTC	Erythromycin	15080-01	1
Euro OTC	Erythromycin	14044205	5
Euro OTC	Erythromycin	150803401	3
Euro OTC	Erythromycin	14090-9-01	1
Euro OTC	Erythromycin	140901902	4
Euro OTC	Erythromycin	1409019-02	55
Euro OTC	Erythromycin	140919-02	1
Euro OTC	Erythromycin	130801602	4
Euro OTC	Erythromycin	181119-01	1
Euro OTC	Erythromycin	1607038-2	1
Euro OTC	Erythromycin	170503803	1
Euro OTC	Erythromycin	17088005	4
Euro OTC	Erythromycin	700770	1
Euro OTC	Erythromycin	1106036-04	2
Euro OTC	Erythromycin	1106036-05	10
Euro OTC	Erythromycin	1705038-03	18
Euro OTC	Erythromycin	95(1705038-02	1
Euro OTC	Erythromycin	1205035-02	19
Euro OTC	Erythromycin	EBS/1205-M-040	1
Euro OTC	Erythromycin	4112Q-01362	1
Euro OTC	Erythromycin	1013u-01362	1
Euro OTC	Erythromycin	17088210	13
Euro OTC	Erythromycin	1811019-01	18
Euro OTC	Erythromycin	170503801	1
Euro OTC	Erythromycin	96(1811019-01)	1
Euro OTC	Erythromycin	1205036-01	11
Euro OTC	Erythromycin	3613A1362	1
Euro OTC	Erythromycin	1205036-02	11
Euro OTC	Erythromycin	1705038-02	10
Euro OTC	Erythromycin	120503602	2
Euro OTC	Erythromycin	1508034	1
Euro OTC	Erythromycin	130801601	6
Euro OTC	Erythromycin	180201E/1705038-01	1
Euro OTC	Erythromycin	EBS/1701-BM-001	3
Euro OTC	Erythromycin	1705038-01	77
Euro OTC	Erythromycin	1205036/02	1
Euro OTC	Erythromycin	160703802	2
Euro OTC	Erythromycin	1308016-02	34
Euro OTC	Erythromycin	1607038/01	1
Euro OTC	Erythromycin	EBS/1307-BM-001	1
Euro OTC	Erythromycin	12A02-N06	1
Euro OTC	Erythromycin	170538-03	1
Euro OTC	Erythromycin	1508034-01	37
Euro OTC	Escin soluble	1612012-01	2
Euro OTC	Escin soluble	1310020-01	1
Euro OTC	Escin soluble	1108038	1
Euro OTC	Estradiol benzoate	130100601	1
Euro OTC	Estradiol benzoate	1412033-02	1
Euro OTC	Estradiol benzoate	1807161	1
Euro OTC	Estradiol benzoate	1611004-028638831	1
Euro OTC	Estradiol benzoate	1611004-02	1
Euro OTC	Estradiol benzoate	1201019-02	12
Euro OTC	Estradiol benzoate	140101-01	1
Euro OTC	Estradiol benzoate	121019-01	1
Euro OTC	Estradiol benzoate	1412033-01	8
Euro OTC	Estradiol benzoate	84303830	2
Euro OTC	Estradiol benzoate	83101340	3
Euro OTC	Estradiol benzoate	1301006-01	3
Euro OTC	Estradiol benzoate	161100402	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Estradiol benzoate	1611004-01	8
Euro OTC	Estradiol benzoate	8011422	1
Euro OTC	Estradiol benzoate	14010001-01	1
Euro OTC	Estradiol benzoate	1401001-01	11
Euro OTC	Estradiol benzoate	1301006-02	11
Euro OTC	Estradiol benzoate	1412033-01	2
Euro OTC	Estradiol benzoate	1007037-02	1
Euro OTC	Estradiol benzoate	150700901	2
Euro OTC	Estradiol benzoate	1507009-01	15
Euro OTC	Estradiol benzoate	1507009-01	1
Euro OTC	Estradiol benzoate	1201019-01	7
Euro OTC	Estradiol hemihydrate	1808032-01	19
Euro OTC	Estradiol hemihydrate	17058205	1
Euro OTC	Estradiol hemihydrate	1601017-01	6
Euro OTC	Estradiol hemihydrate	1607003-01	8
Euro OTC	Estradiol hemihydrate	1404030-01	1
Euro OTC	Estradiol hemihydrate	1504034-01	1
Euro OTC	Estradiol hemihydrate	15D01-B02	1
Euro OTC	Estradiol hemihydrate	1607003-0101374814	1
Euro OTC	Estradiol hemihydrate	1308003-01	2
Euro OTC	Estradiol hemihydrate	1612008-01	26
Euro OTC	Estradiol valerate	15C10-b05	1
Euro OTC	Estriol	1403011-01	1
Euro OTC	Estriol	1403011-01	1
Euro OTC	Estriol	1606044-01	12
Euro OTC	Estriol	1703044-01	1
Euro OTC	Estriol	1610035-01	5
Euro OTC	Estriol	1803011-01	2
Euro OTC	Estriol	15A15-B07-314241	3
Euro OTC	Estriol	1409018-01	1
Euro OTC	Estriol	1201020-01	1
Euro OTC	Estriol	170344-01	1
Euro OTC	Estriol	1310002-01	1
Euro OTC	Estriol	15028904	9
Euro OTC	Estriol	1710030-01	5
Euro OTC	Ethacridine lactate m...	2905007-01	1
Euro OTC	Ethacridine lactate m...	1412034-01	5
Euro OTC	Ethacridine lactate m...	1412034-03	5
Euro OTC	Ethacridine lactate m...	1206016-03	6
Euro OTC	Ethacridine lactate m...	2905007v-03	1
Euro OTC	Ethacridine lactate m...	2905007V-03	1
Euro OTC	Ethacridine lactate m...	1206016-03	1
Euro OTC	Ethacridine lactate m...	1412034-04	1
Euro OTC	Ethacridine lactate m...	1206016-02	1
Euro OTC	Ethacridine lactate m...	1412034-01	1
Euro OTC	Ethacridine lactate m...	1206016-02	3
Euro OTC	Ethacridine lactate m...	1412034-02	2
Euro OTC	Ethyl hydrogen fumara...	1510005-01	12
Euro OTC	Finasteride (form III)	1403016-02	1
Euro OTC	Finasteride (form III)	C1012/116010	1
Euro OTC	Finasteride (form III)	1704004-01	5
Euro OTC	Finasteride (form III)	1605020-01	1
Euro OTC	Finasteride (form III)	1610012-01	2
Euro OTC	Finasteride (form III)	1411008-01	1
Euro OTC	Fluocinololacetone, ...	2804013E-01	1
Euro OTC	Fluocinololacetone, ...	1408022-02	2
Euro OTC	Folic acid hydrate	1604036-01	1
Euro OTC	Folic acid hydrate	1611009-01	4
Euro OTC	Folic acid hydrate	1608030-01	1
Euro OTC	Folic acid hydrate	1611009-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Folic acid hydrate	1401005-01	1
Euro OTC	Folic acid hydrate	1307041-01	1
Euro OTC	Frankincense extract ...	1104005V-04	2
Euro OTC	Fructose	1701049-01	12
Euro OTC	Fructose	1209010-01	1
Euro OTC	Fructose	1212015-01	3
Euro OTC	Fructose	1201017-01	2
Euro OTC	Fructose	1307022-01	5
Euro OTC	Fructose	140100401	1
Euro OTC	Fructose	121201-01pzn	1
Euro OTC	Fructose	1401004-01	5
Euro OTC	Fructose	AD16C91215	1
Euro OTC	Fructose	1710034-01	16
Euro OTC	Fructose	1803040-01	7
Euro OTC	Fructose	180304001	1
Euro OTC	Fructose	141000901	2
Euro OTC	Fructose	AD15190330	1
Euro OTC	Fructose	1803040-1	1
Euro OTC	Fructose	1811018-01	1
Euro OTC	Fructose	1705002-01	7
Euro OTC	Fructose	1803040-012	1
Euro OTC	Fructose	1512009-01	7
Euro OTC	Fructose	1610033-01	11
Euro OTC	Fructose	1410009-01	4
Euro OTC	Gelatine, white	1106021-01	3
Euro OTC	Gentamicin sulfate	FG1501009	1
Euro OTC	Gentamicin sulfate	1305047-03	8
Euro OTC	Gentamicin sulfate	1610036-02	10
Euro OTC	Gentamicin sulfate	180208G/1711025-01	1
Euro OTC	Gentamicin sulfate	1407024-04	5
Euro OTC	Gentamicin sulfate	140702404	4
Euro OTC	Gentamicin sulfate	1711025-01	19
Euro OTC	Gentamicin sulfate	1711025-02	6
Euro OTC	Gentamicin sulfate	80616-2	1
Euro OTC	Gentamicin sulfate	1305047-03	1
Euro OTC	Gentamicin sulfate	1506016-02	27
Euro OTC	Gentamicin sulfate	FG11112482	1
Euro OTC	Gentamicin sulfate	1205021-01	12
Euro OTC	Gentamicin sulfate	150601602	2
Euro OTC	Gentamicin sulfate	4318A-01582	1
Euro OTC	Gentamicin sulfate	4988304	1
Euro OTC	Gentamicin sulfate	1106016-01	5
Euro OTC	Gentamicin sulfate	1604054-03	1
Euro OTC	Gentamicin sulfate	1106016-03	5
Euro OTC	Gentamicin sulfate	1106016-02	3
Euro OTC	Gentamicin sulfate	150601601	1
Euro OTC	Gentamicin sulfate	2911034-02	1
Euro OTC	Gentamicin sulfate	1604054-01	22
Euro OTC	Gentamicin sulfate	17358609	14
Euro OTC	Gentamicin sulfate	1407024-02	1
Euro OTC	Gentamicin sulfate	1610036-015339	1
Euro OTC	Gentamicin sulfate	1407024-03	2
Euro OTC	Gentamicin sulfate	140702403	1
Euro OTC	Gentamicin sulfate	160405401	1
Euro OTC	Gentamicin sulfate	1205021-01	1
Euro OTC	Gentamicin sulfate	130504703	1
Euro OTC	Gentamicin sulfate	1407024-01	5
Euro OTC	Gentamicin sulfate	130504701	1
Euro OTC	Gentamicin sulfate	1205021-02	3
Euro OTC	Gentamicin sulfate	130504703	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Gentamicin sulfate	1205021-02	15
Euro OTC	Gentamicin sulfate	1407024	1
Euro OTC	Gentamicin sulfate	1205021-82	1
Euro OTC	Gentamicin sulfate	1305047-01	3
Euro OTC	Gentamicin sulfate	1106016-03	2
Euro OTC	Gentamicin sulfate	1205021-02	2
Euro OTC	Gentamicin sulfate	1205021-01	1
Euro OTC	Gentamicin sulfate	1305047-03	10
Euro OTC	Ginseng (H1018)	15001047/0	1
Euro OTC	Glucose	18012905	6
Euro OTC	Glucose	18012922	1
Euro OTC	Glucose monohydrate	1410019-01	7
Euro OTC	Glucose monohydrate	140101902	2
Euro OTC	Glucose monohydrate	141001901	1
Euro OTC	Glucose monohydrate	31102013B	1
Euro OTC	Glucose monohydrate	1114U2209	1
Euro OTC	Glucose monohydrate	1305030-02	1
Euro OTC	Glucose monohydrate	1401019-01	19
Euro OTC	Glucose monohydrate	140101901	1
Euro OTC	Glucose monohydrate	2013U-02209	1
Euro OTC	Glucose monohydrate	1405004-02	5
Euro OTC	Glucose monohydrate	1103027-01	1
Euro OTC	Glucose monohydrate	1305030-01	4
Euro OTC	Glucose monohydrate	130503001	1
Euro OTC	Glucose monohydrate	E267X	1
Euro OTC	Glucose monohydrate	1210034-02	11
Euro OTC	Glucose monohydrate	2023012-1	1
Euro OTC	Glucose monohydrate	2313E-02209	1
Euro OTC	Glucose monohydrate	1406212209	1
Euro OTC	Glucose monohydrate	1405004-01	16
Euro OTC	Glucose monohydrate	1204044-02	4
Euro OTC	Glucose monohydrate	140004-01	1
Euro OTC	Glucose monohydrate	14005004-02	1
Euro OTC	Glucose monohydrate	1210034-01	9
Euro OTC	Glucose monohydrate	1807005-02	2
Euro OTC	Glucose monohydrate	190402901	1
Euro OTC	Glucose monohydrate	1701002-01	22
Euro OTC	Glucose monohydrate	1607024	1
Euro OTC	Glucose monohydrate	1504028-01	15
Euro OTC	Glucose monohydrate	1701002-02	6
Euro OTC	Glucose monohydrate	1707004-001	2
Euro OTC	Glucose monohydrate	1707004	1
Euro OTC	Glucose monohydrate	160704-01	1
Euro OTC	Glucose monohydrate	170704401	1
Euro OTC	Glucose monohydrate	180100901	2
Euro OTC	Glucose monohydrate	23051802	1
Euro OTC	Glucose monohydrate	1801009_01	1
Euro OTC	Glucose monohydrate	E617K	1
Euro OTC	Glucose monohydrate	180100-01	1
Euro OTC	Glucose monohydrate	E255L	1
Euro OTC	Glucose monohydrate	1801009-01	39
Euro OTC	Glucose monohydrate	16040060-01	1
Euro OTC	Glucose monohydrate	1504028-02	10
Euro OTC	Glucose monohydrate	1210034-03	1
Euro OTC	Glucose monohydrate	1604006-02	30
Euro OTC	Glucose monohydrate	160400602	2
Euro OTC	Glucose monohydrate	1607024-01	18
Euro OTC	Glucose monohydrate	20160405-04	1
Euro OTC	Glucose monohydrate	1509006-1	1
Euro OTC	Glucose monohydrate	1607024-02	28

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Glucose monohydrate	160101601	2
Euro OTC	Glucose monohydrate	1601016-02	5
Euro OTC	Glucose monohydrate	1601016-01	20
Euro OTC	Glucose monohydrate	1305030-02	6
Euro OTC	Glucose monohydrate	1509006-02	4
Euro OTC	Glucose monohydrate	150402801	1
Euro OTC	Glucose monohydrate	150900601	1
Euro OTC	Glucose monohydrate	15040281	1
Euro OTC	Glucose monohydrate	1509006-01	19
Euro OTC	Glucose monohydrate	1410019-02	5
Euro OTC	Glucose monohydrate	96-01	1
Euro OTC	Glucose monohydrate	170700401	1
Euro OTC	Glucose monohydrate	187005-001	1
Euro OTC	Glucose monohydrate	180100941	1
Euro OTC	Glucose monohydrate	1807005-01	15
Euro OTC	Glucose monohydrate	1811020-01	8
Euro OTC	Glucose monohydrate	1904029-01	1
Euro OTC	Glucose monohydrate	18070521	1
Euro OTC	Glucose monohydrate	18070501	1
Euro OTC	Glucose monohydrate	180700501	1
Euro OTC	Glutamic acid	28071505	1
Euro OTC	Glutamine	1706010-01	1
Euro OTC	Glycine	1801005-01	1
Euro OTC	Glycine (H1021)	16000239/0	2
Euro OTC	Griseofulvin	14E16-B01-300367	1
Euro OTC	Guaifenesin	1706006-01	2
Euro OTC	Guaifenesin	1804004-01	1
Euro OTC	Hamamelis leaf, cut	12164805	1
Euro OTC	Hard fat in pastilles	17107304	3
Euro OTC	Hard fat W35	17K09B03343780	1
Euro OTC	Heparin sodium	1605022-02	1
Euro OTC	Histamine dihydrochlo...	1506030-01	3
Euro OTC	Histamine dihydrochlo...	1111004-01	4
Euro OTC	Histamine dihydrochlo...	111004-02	1
Euro OTC	Histamine dihydrochlo...	1305006-01	1
Euro OTC	Histamine dihydrochlo...	3031503	1
Euro OTC	Histamine dihydrochlo...	20021412	1
Euro OTC	Histamine dihydrochlo...	1111004-02	2
Euro OTC	Histamine dihydrochlo...	1104004-03	1
Euro OTC	Histamine dihydrochlo...	1411012-02	3
Euro OTC	Hyaluronic acid sodiu...	L1512025-01	2
Euro OTC	Hydrocortisone	11611016-01	1
Euro OTC	Hydrocortisone	181102182	1
Euro OTC	Hydrocortisone	1501034	1
Euro OTC	Hydrocortisone	17351202	8
Euro OTC	Hydrocortisone	161101601	2
Euro OTC	Hydrocortisone	1802009-01	23
Euro OTC	Hydrocortisone	161016-01	1
Euro OTC	Hydrocortisone	150902901	1
Euro OTC	Hydrocortisone	16272311	6
Euro OTC	Hydrocortisone	180200901	1
Euro OTC	Hydrocortisone	1611016-01	58
Euro OTC	Hydrocortisone	1405001	1
Euro OTC	Hydrocortisone	715Q-01719	1
Euro OTC	Hydrocortisone	150902	1
Euro OTC	Hydrocortisone	1509029	1
Euro OTC	Hydrocortisone	1505010601	1
Euro OTC	Hydrocortisone	150929-01	1
Euro OTC	Hydrocortisone	1509029-01	42
Euro OTC	Hydrocortisone	157015	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Hydrocortisone	150501601	1
Euro OTC	Hydrocortisone	150516-01	1
Euro OTC	Hydrocortisone	1505016-01	28
Euro OTC	Hydrocortisone	3817A-01719	1
Euro OTC	Hydrocortisone	1107005-01	1
Euro OTC	Hydrocortisone	1107005-02	5
Euro OTC	Hydrocortisone	1611016	2
Euro OTC	Hydrocortisone	4112A-01719	1
Euro OTC	Hydrocortisone	1206018-01	21
Euro OTC	Hydrocortisone	1204307-01	1
Euro OTC	Hydrocortisone	1204037-01	23
Euro OTC	Hydrocortisone	210818A	1
Euro OTC	Hydrocortisone	1408023-01	5
Euro OTC	Hydrocortisone	1107005-03	7
Euro OTC	Hydrocortisone	82(1802009-01)	1
Euro OTC	Hydrocortisone	81(1802009-01)	1
Euro OTC	Hydrocortisone	80(1802009-01)	1
Euro OTC	Hydrocortisone	79(1802009-01)	1
Euro OTC	Hydrocortisone	1802008-01	1
Euro OTC	Hydrocortisone acetate	1706021-01	54
Euro OTC	Hydrocortisone acetate	HCAZ11004K	1
Euro OTC	Hydrocortisone acetate	170602101	2
Euro OTC	Hydrocortisone acetate	Euro0TC	2
Euro OTC	Hydrocortisone acetate	170602-01	1
Euro OTC	Hydrocortisone acetate	1508029-02	1
Euro OTC	Hydrocortisone acetate	17086204	17
Euro OTC	Hydrocortisone acetate	120503401	3
Euro OTC	Hydrocortisone acetate	170830HA/1706021-01	1
Euro OTC	Hydrocortisone acetate	1112031-01	13
Euro OTC	Hydrocortisone acetate	2604052-01	1
Euro OTC	Hydrocortisone acetate	150802901	1
Euro OTC	Hydrocortisone acetate	1107004-03	4
Euro OTC	Hydrocortisone acetate	160405201	5
Euro OTC	Hydrocortisone acetate	2715Q-01718	1
Euro OTC	Hydrocortisone acetate	1508029	2
Euro OTC	Hydrocortisone acetate	911I-01718	1
Euro OTC	Hydrocortisone acetate	1706021	1
Euro OTC	Hydrocortisone acetate	18050101	1
Euro OTC	Hydrocortisone acetate	17066021-01	1
Euro OTC	Hydrocortisone acetate	1409005-01	17
Euro OTC	Hydrocortisone acetate	1409005-02	4
Euro OTC	Hydrocortisone acetate	150103501	3
Euro OTC	Hydrocortisone acetate	1404017-03	1
Euro OTC	Hydrocortisone acetate	1505008-01	46
Euro OTC	Hydrocortisone acetate	140900502	1
Euro OTC	Hydrocortisone acetate	1502000801	1
Euro OTC	Hydrocortisone acetate	1508029-01	32
Euro OTC	Hydrocortisone acetate	L68934	1
Euro OTC	Hydrocortisone acetate	1404017-02	6
Euro OTC	Hydrocortisone acetate	150500801	3
Euro OTC	Hydrocortisone acetate	1505008	1
Euro OTC	Hydrocortisone acetate	1604052-01	72
Euro OTC	Hydrocortisone acetate	508029-01	2
Euro OTC	Hydrocortisone acetate	1404017-01	12
Euro OTC	Hydrocortisone acetate	12090003-01	3
Euro OTC	Hydrocortisone acetate	120900301	2
Euro OTC	Hydrocortisone acetate	1209003-01	31
Euro OTC	Hydrocortisone acetate	180501001	3
Euro OTC	Hydrocortisone acetate	1805010-01	9
Euro OTC	Hydrocortisone acetate	140900501	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Hydrocortisone acetate	1501035-01	14
Euro OTC	Hydroxyethylcellulose...	1108027-01	2
Euro OTC	Hydroxyethylcellulose...	N1535	3
Euro OTC	Hydroxyethylcellulose...	1610036-01	1
Euro OTC	Hydroxyethylcellulose...	1406017	1
Euro OTC	Hydroxyethylcellulose...	1610039-03	9
Euro OTC	Hydroxyethylcellulose...	1500418-02	1
Euro OTC	Hydroxyethylcellulose...	161003902	1
Euro OTC	Hydroxyethylcellulose...	130902003	1
Euro OTC	Hydroxyethylcellulose...	1504019-02	1
Euro OTC	Hydroxyethylcellulose...	3081603	1
Euro OTC	Hydroxyethylcellulose...	914E-03301	1
Euro OTC	Hydroxyethylcellulose...	1406017-03	2
Euro OTC	Hydroxyethylcellulose...	13099020-01	1
Euro OTC	Hydroxyethylcellulose...	160039-01	1
Euro OTC	Hydroxyethylcellulose...	80616-3	1
Euro OTC	Hydroxyethylcellulose...	1711012-01	56
Euro OTC	Hydroxyethylcellulose...	1309020-03	10
Euro OTC	Hydroxyethylcellulose...	180201HC/1610039-01	1
Euro OTC	Hydroxyethylcellulose...	1606030-03	7
Euro OTC	Hydroxyethylcellulose...	1504018-03	2
Euro OTC	Hydroxyethylcellulose...	160603001	3
Euro OTC	Hydroxyethylcellulose...	17319626	2
Euro OTC	Hydroxyethylcellulose...	1907013-01	1
Euro OTC	Hydroxyethylcellulose...	1406017-02	28
Euro OTC	Hydroxyethylcellulose...	1309020-01	1
Euro OTC	Hydroxyethylcellulose...	1496017-02	1
Euro OTC	Hydroxyethylcellulose...	1406017-02	15
Euro OTC	Hydroxyethylcellulose...	13072015A	1
Euro OTC	Hydroxyethylcellulose...	140601702	2
Euro OTC	Hydroxyethylcellulose...	160603	1
Euro OTC	Hydroxyethylcellulose...	1309020-02	6
Euro OTC	Hydroxyethylcellulose...	1105054-01	1
Euro OTC	Hydroxyethylcellulose...	1504018-01	15
Euro OTC	Hydroxyethylcellulose...	1504018-02	30
Euro OTC	Hydroxyethylcellulose...	1406017-01	8
Euro OTC	Hydroxyethylcellulose...	41(1907013-01)	1
Euro OTC	Hydroxyethylcellulose...	16100603	2
Euro OTC	Hydroxyethylcellulose...	1606030-01	2
Euro OTC	Hydroxyethylcellulose...	154018-02	1
Euro OTC	Hydroxyethylcellulose...	1504018-02	2
Euro OTC	Hydroxyethylcellulose...	1504018-01	1
Euro OTC	Hydroxyethylcellulose...	12175708	1
Euro OTC	Hydroxyethylcellulose...	140601703	1
Euro OTC	Hydroxyethylcellulose...	1110024-04-336431909	1
Euro OTC	Hydroxyethylcellulose...	P0745	3
Euro OTC	Hydroxyethylcellulose...	1711012	1
Euro OTC	Hydroxyethylcellulose...	110505401	1
Euro OTC	Hydroxyethylcellulose...	181102185	1
Euro OTC	Hydroxyethylcellulose...	16603001	1
Euro OTC	Hydroxyethylcellulose...	150401801	1
Euro OTC	Hydroxyethylcellulose...	1504018-03	21
Euro OTC	Hydroxyethylcellulose...	130902001	2
Euro OTC	Hydroxyethylcellulose...	1406212528	1
Euro OTC	Hydroxyethylcellulose...	17042014A	1
Euro OTC	Hydroxyethylcellulose...	130400801	4
Euro OTC	Hydroxyethylcellulose...	40(1711012-01)	1
Euro OTC	Hydroxyethylcellulose...	16060302-02	1
Euro OTC	Hydroxyethylcellulose...	161003901	4
Euro OTC	Hydroxyethylcellulose...	1606030-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Hydroxyethylcellulose...	1108026-01	1
Euro OTC	Hydroxyethylcellulose...	K1371	2
Euro OTC	Hydroxyethylcellulose...	1309020-01	25
Euro OTC	Hydroxyethylcellulose...	111024-02	1
Euro OTC	Hydroxyethylcellulose...	13040008-02	1
Euro OTC	Hydroxyethylcellulose...	1110024-02	7
Euro OTC	Hydroxyethylcellulose...	1613M-03301	1
Euro OTC	Hydroxyethylcellulose...	130400802	2
Euro OTC	Hydroxyethylcellulose...	1108026-01	1
Euro OTC	Hydroxyethylcellulose...	1007018-01	1
Euro OTC	Hydroxyethylcellulose...	1610039-01	53
Euro OTC	Hydroxyethylcellulose...	1110024-04	22
Euro OTC	Hydroxyethylcellulose...	1304008-02	23
Euro OTC	Hydroxyethylcellulose...	1610039-01	1
Euro OTC	Hydroxyethylcellulose...	1406017-01	1
Euro OTC	Hydroxyethylcellulose...	11110024-03	1
Euro OTC	Hydroxyethylcellulose...	1504018-03	1
Euro OTC	Hydroxyethylcellulose...	16012909	6
Euro OTC	Hydroxyethylcellulose...	314Q3301	1
Euro OTC	Hydroxyethylcellulose...	146017-02	1
Euro OTC	Ibuprofen	9061807	1
Euro OTC	Ibuprofen	1410025-02	7
Euro OTC	Ibuprofen	1410025-01	6
Euro OTC	Ibuprofen	1211011-01	1
Euro OTC	Indometacin	1801040-01	6
Euro OTC	Indometacin	1612017-02	3
Euro OTC	Indometacin	1612017-01	2
Euro OTC	Indometacin	180104002	1
Euro OTC	Indometacin	1801040-02	2
Euro OTC	Indometacin	201606272030	1
Euro OTC	Indometacin	180104001	1
Euro OTC	Indometacin	1512031-01	9
Euro OTC	Indometacin	1205002-01	2
Euro OTC	Inositol	1609025-01	1
Euro OTC	Inositol	130400901	1
Euro OTC	Inositol	1304009-01	1
Euro OTC	Isosorbide dinitrate ...	15I10-B01325474	1
Euro OTC	Isosorbide dinitrate ...	PN:1370-240409	1
Euro OTC	Lactose monohydrate	160602501	1
Euro OTC	Lactose monohydrate	1606025-01	27
Euro OTC	Lactose monohydrate	1301028-03	2
Euro OTC	Lactose monohydrate	15090-30-02	1
Euro OTC	Lactose monohydrate	150101701	1
Euro OTC	Lactose monohydrate	1401003	1
Euro OTC	Lactose monohydrate	15010	1
Euro OTC	Lactose monohydrate	130102801	1
Euro OTC	Lactose monohydrate	1301028-02	17
Euro OTC	Lactose monohydrate	31102013A	1
Euro OTC	Lactose monohydrate	L1237A4172	2
Euro OTC	Lactose monohydrate	1301028-01	9
Euro OTC	Lactose monohydrate	1401003-01	14
Euro OTC	Lactose monohydrate	110602202	1
Euro OTC	Lactose monohydrate	1112020-02	9
Euro OTC	Lactose monohydrate	140-1003-01	1
Euro OTC	Lactose monohydrate	18082014A	1
Euro OTC	Lactose monohydrate	1401003-02	7
Euro OTC	Lactose monohydrate	140100302	1
Euro OTC	Lactose monohydrate	5123	1
Euro OTC	Lactose monohydrate	1501017-02	1
Euro OTC	Lactose monohydrate	1509030-01	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Lactose monohydrate	3215A-02210	1
Euro OTC	Lactose monohydrate	1013Q-02210	1
Euro OTC	Lactose monohydrate	1509030-02	9
Euro OTC	Lactose monohydrate	1802002-01	9
Euro OTC	Lactose monohydrate	L101854417A545	1
Euro OTC	Lactose monohydrate	UP250171017	1
Euro OTC	Lactose monohydrate	1812008-01	1
Euro OTC	Lactose monohydrate	180200201	1
Euro OTC	Lactose monohydrate	L101852618	1
Euro OTC	Lactose monohydrate	1606025-02	3
Euro OTC	Levocarnitine	L1801028-02	1
Euro OTC	Levocarnitine	1303027	1
Euro OTC	Levocarnitine	L1705006-01	1
Euro OTC	Levocarnitine	L1608038-01	2
Euro OTC	Levocarnitine	1705006-01	1
Euro OTC	Levothyroxin sodium	17J10-B05-345570	1
Euro OTC	Lidocaine	161202102	2
Euro OTC	Lidocaine	1603027-01	2
Euro OTC	Lidocaine	1310023-02	3
Euro OTC	Lidocaine	1612021-01	1
Euro OTC	Lidocaine hydrochlori...	1202035/02	1
Euro OTC	Lidocaine hydrochlori...	20131010-4	1
Euro OTC	Lidocaine hydrochlori...	1606021-01	1
Euro OTC	Lidocaine hydrochlori...	16006024-02	1
Euro OTC	Lidocaine hydrochlori...	1606024-01	1
Euro OTC	Lidocaine hydrochlori...	1202035-01	28
Euro OTC	Lidocaine hydrochlori...	140415A	1
Euro OTC	Lidocaine hydrochlori...	2810029-01	1
Euro OTC	Lidocaine hydrochlori...	160602401	4
Euro OTC	Lidocaine hydrochlori...	1202035-1	1
Euro OTC	Lidocaine hydrochlori...	2213U01775	1
Euro OTC	Lidocaine hydrochlori...	1808031-01	31
Euro OTC	Lidocaine hydrochlori...	23081601	1
Euro OTC	Lidocaine hydrochlori...	1505023-03	19
Euro OTC	Lidocaine hydrochlori...	1411010-01	1
Euro OTC	Lidocaine hydrochlori...	1505023-01	63
Euro OTC	Lidocaine hydrochlori...	141101001	1
Euro OTC	Lidocaine hydrochlori...	1411010-2	1
Euro OTC	Lidocaine hydrochlori...	1108017-01	1
Euro OTC	Lidocaine hydrochlori...	1202034-01	2
Euro OTC	Lidocaine hydrochlori...	1411010-02	12
Euro OTC	Lidocaine hydrochlori...	1411010-02Mg50g	1
Euro OTC	Lidocaine hydrochlori...	141101002	3
Euro OTC	Lidocaine hydrochlori...	1606024-01	78
Euro OTC	Lidocaine hydrochlori...	1411010-02	13
Euro OTC	Lidocaine hydrochlori...	150502301	2
Euro OTC	Lidocaine hydrochlori...	14110101-01	1
Euro OTC	Lidocaine hydrochlori...	141101001	3
Euro OTC	Lidocaine hydrochlori...	1505023- ~~8~~ 1	1
Euro OTC	Lidocaine hydrochlori...	1401010-01	1
Euro OTC	Lidocaine hydrochlori...	141101-01	1
Euro OTC	Lidocaine hydrochlori...	141101002	1
Euro OTC	Lidocaine hydrochlori...	141101	1
Euro OTC	Lidocaine hydrochlori...	1505023-02	21
Euro OTC	Lidocaine hydrochlori...	1411010-01Mg25g	1
Euro OTC	Lidocaine hydrochlori...	1066024-01	1
Euro OTC	Lidocaine hydrochlori...	1403025-03	1
Euro OTC	Lidocaine hydrochlori...	160602402	1
Euro OTC	Lidocaine hydrochlori...	12020351612015-01	1
Euro OTC	Lidocaine hydrochlori...	1606024-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Lidocaine hydrochlori...	1612015-02	48
Euro OTC	Lidocaine hydrochlori...	16112015-04	1
Euro OTC	Lidocaine hydrochlori...	1612015	4
Euro OTC	Lidocaine hydrochlori...	180220L/1612015-01	1
Euro OTC	Lidocaine hydrochlori...	1808031	1
Euro OTC	Lidocaine hydrochlori...	ldh/116023	1
Euro OTC	Lidocaine hydrochlori...	1103034-01	2
Euro OTC	Lidocaine hydrochlori...	1808031-02	6
Euro OTC	Lidocaine hydrochlori...	1808031-03	1
Euro OTC	Lidocaine hydrochlori...	1618015-02	1
Euro OTC	Lidocaine hydrochlori...	1612015-03	11
Euro OTC	Lidocaine hydrochlori...	161201504	2
Euro OTC	Lidocaine hydrochlori...	1612015-04	13
Euro OTC	Lidocaine hydrochlori...	1612015-05	1
Euro OTC	Lidocaine hydrochlori...	1612012015-04	1
Euro OTC	Lidocaine hydrochlori...	1202035-03	16
Euro OTC	Lidocaine hydrochlori...	16112015-01	2
Euro OTC	Lidocaine hydrochlori...	171106LHC/1606024-03	1
Euro OTC	Lidocaine hydrochlori...	180803101	2
Euro OTC	Lidocaine hydrochlori...	170612182	1
Euro OTC	Lidocaine hydrochlori...	140302502	2
Euro OTC	Lidocaine hydrochlori...	1505023-01	10
Euro OTC	Lidocaine hydrochlori...	1600024-02	1
Euro OTC	Lidocaine hydrochlori...	1612015-01	41
Euro OTC	Lidocaine hydrochlori...	1411010-01	37
Euro OTC	Lidocaine hydrochlori...	1403025-02Mg:25g	1
Euro OTC	Lidocaine hydrochlori...	13343004Mg:10g	1
Euro OTC	Lidocaine hydrochlori...	LID-20-11-001	2
Euro OTC	Lidocaine hydrochlori...	1403025-02	1
Euro OTC	Lidocaine hydrochlori...	1403025-01	21
Euro OTC	Lidocaine hydrochlori...	1606024-03	13
Euro OTC	Lidocaine hydrochlori...	1004040-01	1
Euro OTC	Lidocaine hydrochlori...	1202035	2
Euro OTC	Lidocaine hydrochlori...	2805016-01	1
Euro OTC	Lidocaine hydrochlori...	1010024-01	2
Euro OTC	Lidocaine hydrochlori...	1108017-02	8
Euro OTC	Lidocaine hydrochlori...	120203502	3
Euro OTC	Lidocaine hydrochlori...	1140917/160602403	1
Euro OTC	L-isoleucine	1708027-01	1
Euro OTC	L-isoleucine	1704022-01	1
Euro OTC	L-leucine	1312032-01	1
Euro OTC	L-leucine	1708010-01	2
Euro OTC	L-phenylalanine	1704003-02	1
Euro OTC	L-phenylalanine	1409013-02	12
Euro OTC	L-proline	L1210046-03	1
Euro OTC	L-threonine	1606019-02	1
Euro OTC	L-tryptophan	150600302	1
Euro OTC	L-tryptophan	1307004-01	1
Euro OTC	L-tryptophan	1307004-01	1
Euro OTC	L-tryptophan	1506003-03	2
Euro OTC	L-tryptophan	1712006-01	1
Euro OTC	L-tryptophan	180600601	1
Euro OTC	Lysine hydrochloride	1805030-01	1
Euro OTC	Lysine hydrochloride	1708009-01	2
Euro OTC	Macrogols 1500	1103017-02	4
Euro OTC	Macrogols 1500	1503026-01	4
Euro OTC	Macrogols 1500	1206005-01	5
Euro OTC	Macrogols 1500	1406004-01	2
Euro OTC	Macrogols 1500	1304023-01	7
Euro OTC	Macrogols 4000	1305020-02	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Macrogols 4000	130502002	1
Euro OTC	Macrogols 4000	1202040-02	3
Euro OTC	Macrogols 4000	1808003-01	1
Euro OTC	Macrogols 4000	1005011-01	1
Euro OTC	Macrogols 4000	1305020-02	2
Euro OTC	Macrogols 4000	140701001	1
Euro OTC	Macrogols 4000	1508031-01	1
Euro OTC	Macrogols 4000	1507004-01	1
Euro OTC	Macrogols 4000	1508031-02	2
Euro OTC	Macrogols 4000	1810003-01	1
Euro OTC	Magnesium carbonate, ...	1211005-01	4
Euro OTC	Magnesium carbonate, ...	1305024-01	1
Euro OTC	Magnesium chloride he...	1212013-01	3
Euro OTC	Magnesium chloride he...	1204017-01	1
Euro OTC	Magnesium chloride he...	1112013-01	2
Euro OTC	Magnesium citrate	1602028-01	1
Euro OTC	Magnesium citrate	1808020-01	1
Euro OTC	Magnesium citrate	1706007-01	2
Euro OTC	Magnesium hydrogen as...	1605031-01	1
Euro OTC	Magnesium hydrogen ci...	1080144	1
Euro OTC	Magnesium hydrogen ci...	1407048-01	4
Euro OTC	Magnesium hydrogen ci...	140704801	1
Euro OTC	Magnesium hydrogen ci...	16000456/0	1
Euro OTC	Magnesium hydrogen ci...	1710001-01	1
Euro OTC	Magnesium hydrogen ci...	1507010-02	2
Euro OTC	Magnesium hydrogen ci...	1802026-01	5
Euro OTC	Magnesium hydrogen ci...	1510033-01	3
Euro OTC	Magnesium hydrogen ci...	151003301	2
Euro OTC	Magnesium hydrogen ci...	1904011-01	1
Euro OTC	Magnesium hydrogen ci...	1085730	2
Euro OTC	Magnesium hydrogen ci...	1501026-01	3
Euro OTC	Magnesium hydrogen ci...	1808006-01	2
Euro OTC	Magnesium hydrogen ci...	1507010-01	12
Euro OTC	Magnesium hydrogen ci...	1407048-02	3
Euro OTC	Magnesium hydrogen ci...	1604002-01	36
Euro OTC	Magnesium hydrogen ci...	120800602	1
Euro OTC	Magnesium hydrogen ci...	1710007-01	8
Euro OTC	Magnesium hydrogen ci...	160100501	2
Euro OTC	Magnesium hydrogen ci...	1604002-012x	1
Euro OTC	Magnesium hydrogen ci...	1110026-01	2
Euro OTC	Magnesium hydrogen ci...	1208006-01	3
Euro OTC	Magnesium hydrogen ci...	1211033-01	6
Euro OTC	Magnesium hydrogen ci...	1604002-02	7
Euro OTC	Magnesium hydrogen ci...	1211033-02	4
Euro OTC	Magnesium hydrogen ci...	1306034-01	3
Euro OTC	Magnesium hydrogen ci...	1306034-02	2
Euro OTC	Magnesium hydrogen ci...	130603402	1
Euro OTC	Magnesium hydrogen ci...	1311003-01	10
Euro OTC	Magnesium hydrogen ci...	2114M1829	1
Euro OTC	Magnesium hydrogen ci...	1311003-02	2
Euro OTC	Magnesium hydrogen ci...	1608004-01	24
Euro OTC	Magnesium hydrogen ci...	1604002.01	1
Euro OTC	Magnesium hydrogen ci...	1404023-01	5
Euro OTC	Magnesium hydrogen ci...	1705010-01	5
Euro OTC	Magnesium hydrogen ci...	1705009-01	5
Euro OTC	Magnesium hydrogen ci...	1701004-01	4
Euro OTC	Magnesium hydrogen ci...	1608004-02	2
Euro OTC	Magnesium hydrogen ci...	1604002-001	1
Euro OTC	Magnesium hydrogen ci...	1811038-01	1
Euro OTC	Magnesium hydrogen ci...	1701003-01	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Magnesium hydrogen ci...	140402301	1
Euro OTC	Magnesium hydrogen ci...	160800401	2
Euro OTC	Magnesium hydrogen ci...		1
Euro OTC	Magnesium orotate dihydro...	1707018-01	3
Euro OTC	Magnesium orotate dihydro...	151102101	1
Euro OTC	Magnesium orotate dihydro...	1504002-01	1
Euro OTC	Magnesium orotate dihydro...	16085017-01	1
Euro OTC	Magnesium orotate dihydro...	1511021-02	6
Euro OTC	Magnesium orotate dihydro...	1608017-01	2
Euro OTC	Magnesium orotate dihydro...	1802018-01	1
Euro OTC	Magnesium orotate dihydro...	130501102	1
Euro OTC	Magnesium orotate dihydro...	1305011-02	5
Euro OTC	Magnesium orotate dihydro...	1305011-01	1
Euro OTC	Magnesium orotate dihydro...	1209029-01	3
Euro OTC	Magnesium orotate dihydro...	1712001-01	2
Euro OTC	Magnesium orotate dihydro...	1407005-01	2
Euro OTC	Magnesium stearate	1912I-01834	1
Euro OTC	Magnesium stearate	914I-01834	1
Euro OTC	Magnesium stearate	315E-01834	1
Euro OTC	Magnesium stearate	1116E-01834	1
Euro OTC	Magnesium sulfate heptahydr...	170610	1
Euro OTC	Magnesium sulfate heptahydr...	1306015-01	1
Euro OTC	Magnesium sulfate heptahydr...	130300801	1
Euro OTC	Magnesium sulfate heptahydr...	1202017-01	6
Euro OTC	Magnesium sulfate heptahydr...	1109023-01	2
Euro OTC	Mannitol	1409016-1	1
Euro OTC	Mannitol	13189216	2
Euro OTC	Mannitol	1511024-01	1
Euro OTC	Mannitol	180219M/1708005-01	1
Euro OTC	Mannitol	150801	1
Euro OTC	Mannitol	1708005-02	3
Euro OTC	Mannitol	130504301	2
Euro OTC	Mannitol	1802004-01	10
Euro OTC	Mannitol	170304302	1
Euro OTC	Mannitol	E976E	1
Euro OTC	Mannitol	1006019-01	1
Euro OTC	Mannitol	1305043-02	17
Euro OTC	Mannitol	1310004-01	27
Euro OTC	Mannitol	151105-01	1
Euro OTC	Mannitol	17/14	1
Euro OTC	Mannitol	1610027-02	1
Euro OTC	Mannitol	151102501	2
Euro OTC	Mannitol	180200401	1
Euro OTC	Mannitol	1511025-01	39
Euro OTC	Mannitol	158012-01	1
Euro OTC	Mannitol	1502022/01	1
Euro OTC	Mannitol	1806024-01	1
Euro OTC	Mannitol	1806024/01	1
Euro OTC	Mannitol	1704010-01	1
Euro OTC	Mannitol	e118t	1
Euro OTC	Mannitol	150801201	1
Euro OTC	Mannitol	E118T	1
Euro OTC	Mannitol	1902011-01	1
Euro OTC	Mannitol	1809008-02	1
Euro OTC	Mannitol	180900802	1
Euro OTC	Mannitol	140403402	1
Euro OTC	Mannitol	1508012-01	39
Euro OTC	Mannitol	1305043	1
Euro OTC	Mannitol	16006031-01	1
Euro OTC	Mannitol	1606031-02	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Mannitol	161002701	3
Euro OTC	Mannitol	160027-01	1
Euro OTC	Mannitol	16H29-B02	1
Euro OTC	Mannitol	1703043-01	14
Euro OTC	Mannitol	160603101	1
Euro OTC	Mannitol	16040008-01	1
Euro OTC	Mannitol	160400	2
Euro OTC	Mannitol	1104028-03	3
Euro OTC	Mannitol	1703049-01	1
Euro OTC	Mannitol	1110009-01	10
Euro OTC	Mannitol	1104028-01	1
Euro OTC	Mannitol	1205016-01	15
Euro OTC	Mannitol	1703043-02	5
Euro OTC	Mannitol	161129M/1604008-01	1
Euro OTC	Mannitol	160601-01	1
Euro OTC	Mannitol	12105433	9
Euro OTC	Mannitol	E671Z	1
Euro OTC	Mannitol	1604002-01	1
Euro OTC	Mannitol	1205016	1
Euro OTC	Mannitol	150916M	1
Euro OTC	Mannitol	1211034-01	12
Euro OTC	Mannitol	160408-01	1
Euro OTC	Mannitol	1604008	1
Euro OTC	Mannitol	1708005-01	28
Euro OTC	Mannitol	111100009-02	1
Euro OTC	Mannitol	1305043-01	15
Euro OTC	Mannitol	170304301	2
Euro OTC	Mannitol	1606031-01	18
Euro OTC	Mannitol	150202202	1
Euro OTC	Mannitol	1610027-01	31
Euro OTC	Mannitol	1502022-01	20
Euro OTC	Mannitol	1409016-02	7
Euro OTC	Mannitol	1404036-02	1
Euro OTC	Mannitol	1404034-02	17
Euro OTC	Mannitol	10112014E	1
Euro OTC	Mannitol	1409016-01	21
Euro OTC	Mannitol	1404034-01	20
Euro OTC	Mannitol	14003402	1
Euro OTC	Mannitol	T03-07/1409016-01	1
Euro OTC	Mannitol	15002022-01	1
Euro OTC	Mannitol	1310004-02	10
Euro OTC	Mannitol	150202201	2
Euro OTC	Mannitol	140403401	1
Euro OTC	Mannitol	13155710	1
Euro OTC	Marian thistle herb, ...	3701A-111207	1
Euro OTC	Melatonin	1610040-01	2
Euro OTC	Melatonin	1610040-02	6
Euro OTC	Melatonin	1407028-01	1
Euro OTC	Melatonin	120008	1
Euro OTC	Melatonin	180401501	1
Euro OTC	Melatonin	1407028-02	1
Euro OTC	Melatonin	151004101	3
Euro OTC	Melatonin	1510041-01	6
Euro OTC	Melatonin	1407028-02	1
Euro OTC	Melatonin	12001	2
Euro OTC	Melatonin	120400302	1
Euro OTC	Melatonin	1804015-01	3
Euro OTC	Melatonin	1305029-01	1
Euro OTC	Melatonin	1709020-01	12
Euro OTC	Melatonin	1204003-02	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Menthol	1001505	2
Euro OTC	Menthol	1210017-02	11
Euro OTC	Menthol	1210017-01	5
Euro OTC	Menthol	1802003/01	1
Euro OTC	Menthol	1701014-2	1
Euro OTC	Menthol	1310J-01855	1
Euro OTC	Menthol	696773	1
Euro OTC	Menthol	1705022011	1
Euro OTC	Menthol	14030222-02	1
Euro OTC	Menthol	120403602	1
Euro OTC	Menthol	1008029-01	1
Euro OTC	Menthol	1705022	1
Euro OTC	Menthol	1410004-01	18
Euro OTC	Menthol	120-4036-02	2
Euro OTC	Menthol	25021709	1
Euro OTC	Menthol	160801201	2
Euro OTC	Menthol	1701014-01	14
Euro OTC	Menthol	1608012-12	1
Euro OTC	Menthol	20116	1
Euro OTC	Menthol	1701014-02	22
Euro OTC	Menthol	1204036-02	4
Euro OTC	Menthol	170101402	4
Euro OTC	Menthol	1608012-1	1
Euro OTC	Menthol	973192	1
Euro OTC	Menthol	180200301	2
Euro OTC	Menthol	18121718	1
Euro OTC	Menthol	1111020-01	1
Euro OTC	Menthol	1204036-01	11
Euro OTC	Menthol	170101401	1
Euro OTC	Menthol	1204036-1	1
Euro OTC	Menthol	1802003-01	20
Euro OTC	Menthol	1105028-01	1
Euro OTC	Menthol	1305010-02	7
Euro OTC	Menthol	1602024-01	1
Euro OTC	Menthol	1608012-01	19
Euro OTC	Menthol	1603021-01	1
Euro OTC	Menthol	1603024-01	40
Euro OTC	Menthol	150904502	1
Euro OTC	Menthol	1509045-01	4
Euro OTC	Menthol	1509045	2
Euro OTC	Menthol	1502006	1
Euro OTC	Menthol	15020006-02	1
Euro OTC	Menthol	130501002	1
Euro OTC	Menthol	1807044-02	1
Euro OTC	Menthol	1002019-03	1
Euro OTC	Menthol	807862	1
Euro OTC	Menthol	858607	2
Euro OTC	Menthol	1502006-02	17
Euro OTC	Menthol	150200601	1
Euro OTC	Menthol	1502006-01	16
Euro OTC	Menthol	1410004-02	9
Euro OTC	Menthol	180704401	1
Euro OTC	Menthol	141000401	1
Euro OTC	Menthol	140322-01	1
Euro OTC	Menthol	1403022-01	3
Euro OTC	Menthol	6102014B	1
Euro OTC	Menthol	721123	1
Euro OTC	Menthol	140302201	2
Euro OTC	Menthol	PM14	1
Euro OTC	Mesalazine	1602005-02	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Mesalazine	1602005-01	2
Euro OTC	Mesalazine	1206029-04	1
Euro OTC	Methionine	1602012-01	1
Euro OTC	Methoxsalen	1802021-01	6
Euro OTC	Methoxsalen	1608019-01	3
Euro OTC	Methoxsalen	1503009-02	7
Euro OTC	Methoxsalen	1506013-02	7
Euro OTC	Methoxsalen	160801902	1
Euro OTC	Methoxsalen	140900302	2
Euro OTC	Methoxsalen	1503009-01	17
Euro OTC	Methoxsalen	1506013-01	15
Euro OTC	Methoxsalen	120901501	1
Euro OTC	Methoxsalen	1306041-01	5
Euro OTC	Methoxsalen	1203009-01	1
Euro OTC	Methoxsalen	1608019-02	9
Euro OTC	Methoxsalen	1209015-01	7
Euro OTC	Methoxsalen	1104010-01	2
Euro OTC	Methyl-4-hydroxybenzo...	1401027-02	2
Euro OTC	Methyl-4-hydroxybenzo...	1502009-01	8
Euro OTC	Methyl-4-hydroxybenzo...	1009022-01	1
Euro OTC	Methyl-4-hydroxybenzo...	1507042-01	10
Euro OTC	Methyl-4-hydroxybenzo...	2606022	1
Euro OTC	Methyl-4-hydroxybenzo...	172620	1
Euro OTC	Methyl-4-hydroxybenzo...	EA2781	1
Euro OTC	Methyl-4-hydroxybenzo...	1706005-01	1
Euro OTC	Methyl-4-hydroxybenzo...	1758101	1
Euro OTC	Methyl-4-hydroxybenzo...	4314U-01930	1
Euro OTC	Methyl-4-hydroxybenzo...	1507042-02	2
Euro OTC	Methyl-4-hydroxybenzo...	280903701	2
Euro OTC	Methyl-4-hydroxybenzo...	3115M-01930	1
Euro OTC	Methyl-4-hydroxybenzo...	1601037-01	8
Euro OTC	Methyl-4-hydroxybenzo...	1401027-01	1
Euro OTC	Methyl-4-hydroxybenzo...	1807009-01	1
Euro OTC	Methyl-4-hydroxybenzo...	4312I-01930	1
Euro OTC	Methyl-4-hydroxybenzo...	140600501	1
Euro OTC	Methyl-4-hydroxybenzo...	1406005-01	1
Euro OTC	Methylcellulose	16284003	6
Euro OTC	Methylcellulose	10247609	1
Euro OTC	Methylcobalamin	1601043-02	1
Euro OTC	Methylcobalamin	1804006-01	1
Euro OTC	Methylcobalamin	1711031-01	1
Euro OTC	Metronidazole	1109030-01	7
Euro OTC	Metronidazole	1707008-01	27
Euro OTC	Metronidazole	1307042-01	11
Euro OTC	Metronidazole	1508004-01	25
Euro OTC	Metronidazole	170103901	1
Euro OTC	Metronidazole	169028-01	1
Euro OTC	Metronidazole	1707088-01	1
Euro OTC	Metronidazole	1412016-02	21
Euro OTC	Metronidazole	14120601	1
Euro OTC	Metronidazole	13070002-02	1
Euro OTC	Metronidazole	1112028-01	5
Euro OTC	Metronidazole	1407003	1
Euro OTC	Metronidazole	1412016-01	18
Euro OTC	Metronidazole	140601001	3
Euro OTC	Metronidazole	140700302	2
Euro OTC	Metronidazole	1307042-02	7
Euro OTC	Metronidazole	1612013-02	7
Euro OTC	Metronidazole	1407003-02	8
Euro OTC	Metronidazole	130704201	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Metronidazole	1307002-01	14
Euro OTC	Metronidazole	1505024-01	12
Euro OTC	Metronidazole	1407003-01	11
Euro OTC	Metronidazole	1406010-01Mg:50g	1
Euro OTC	Metronidazole	1109030-02	2
Euro OTC	Metronidazole	1212030-02	7
Euro OTC	Metronidazole	1112028-02	8
Euro OTC	Metronidazole	1308704201	1
Euro OTC	Metronidazole	1701039	1
Euro OTC	Metronidazole	150800	1
Euro OTC	Metronidazole	F120017	1
Euro OTC	Metronidazole	176022-01	1
Euro OTC	Metronidazole	170700801	1
Euro OTC	Metronidazole	121001-02	3
Euro OTC	Metronidazole	1508004	1
Euro OTC	Metronidazole	1505024-02	8
Euro OTC	Metronidazole	1707008-02	3
Euro OTC	Metronidazole	1703028-01	1
Euro OTC	Metronidazole	1109030-03	3
Euro OTC	Metronidazole	161201301	2
Euro OTC	Metronidazole	C04-160812	1
Euro OTC	Metronidazole	150804-01	1
Euro OTC	Metronidazole	1701039-01	44
Euro OTC	Metronidazole	15a20-b07-303294	1
Euro OTC	Metronidazole	1508004-02	11
Euro OTC	Metronidazole	1412016	1
Euro OTC	Metronidazole	1612013-01	22
Euro OTC	Metronidazole	1802007-01	5
Euro OTC	Metronidazole	700837	1
Euro OTC	Metronidazole	1901042-01	1
Euro OTC	Metronidazole	1109016-01	1
Euro OTC	Metronidazole	1606028	1
Euro OTC	Metronidazole	1609023-01	1
Euro OTC	Metronidazole	1609028-02	12
Euro OTC	Metronidazole	1802007-02	4
Euro OTC	Metronidazole	150502402	1
Euro OTC	Metronidazole	10800401	1
Euro OTC	Metronidazole	1212001-01	8
Euro OTC	Metronidazole	1212030-01	5
Euro OTC	Metronidazole	1503004-01	1
Euro OTC	Metronidazole	1505024-03	4
Euro OTC	Metronidazole	150502401	1
Euro OTC	Metronidazole	150800401	1
Euro OTC	Metronidazole	20111025	1
Euro OTC	Metronidazole	300915ME	1
Euro OTC	Miconazole nitrate	1611024-015307	1
Euro OTC	Miconazole nitrate	1509014-01	6
Euro OTC	Miconazole nitrate	1811015-01	4
Euro OTC	Miconazole nitrate	1705039-02	2
Euro OTC	Miconazole nitrate	1307025-01	2
Euro OTC	Miconazole nitrate	1611024-015330	1
Euro OTC	Miconazole nitrate	17C01-B02-341814	1
Euro OTC	Miconazole nitrate	1504003-01	2
Euro OTC	Miconazole nitrate	1604004-01	4
Euro OTC	Miconazole nitrate	1705039-015371	1
Euro OTC	Miconazole nitrate	1803028-01	4
Euro OTC	Miconazole nitrate	1705039-01	10
Euro OTC	Miconazole nitrate	1409023-02	4
Euro OTC	Miconazole nitrate	140902302	1
Euro OTC	Miconazole nitrate	1705039-015432	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Miconazole nitrate	1307025-02	3
Euro OTC	Miconazole nitrate	1504903-01	1
Euro OTC	Miconazole nitrate	1611024-02	2
Euro OTC	Miconazole nitrate	1409023-01	1
Euro OTC	Miconazole nitrate	1710020-015457	1
Euro OTC	Miconazole nitrate	17C01-B02-339699	1
Euro OTC	Miconazole nitrate	1607004-01	7
Euro OTC	Miconazole nitrate	1710020-01	8
Euro OTC	Miconazole nitrate	1611024-01	7
Euro OTC	Minoxidil	103004-01	1
Euro OTC	Minoxidil	1504013-02	2
Euro OTC	Minoxidil	1602016-01	1
Euro OTC	Minoxidil	170703M/1602013-02	1
Euro OTC	Minoxidil	12/17	1
Euro OTC	Minoxidil	104029-02	1
Euro OTC	Minoxidil	15081829	1
Euro OTC	Minoxidil	1803004	1
Euro OTC	Minoxidil	1203035-01	12
Euro OTC	Minoxidil	150401-01	1
Euro OTC	Minoxidil	1602013-01	3
Euro OTC	Minoxidil	1710043_01	1
Euro OTC	Minoxidil	1602013-1	1
Euro OTC	Minoxidil	1504013-2	1
Euro OTC	Minoxidil	1607007	1
Euro OTC	Minoxidil	131103801	4
Euro OTC	Minoxidil	160201301	1
Euro OTC	Minoxidil	1919	1
Euro OTC	Minoxidil	1203035-01	1
Euro OTC	Minoxidil	1504013-01	21
Euro OTC	Minoxidil	1203035-03	14
Euro OTC	Minoxidil	1311038-01	38
Euro OTC	Minoxidil	1203035-02	7
Euro OTC	Minoxidil	1607007-01	35
Euro OTC	Minoxidil	2091714	1
Euro OTC	Minoxidil	29071402	1
Euro OTC	Minoxidil	1807018-02	1
Euro OTC	Minoxidil	190100501	1
Euro OTC	Minoxidil	1108039-01	1
Euro OTC	Minoxidil	1411006-01	7
Euro OTC	Minoxidil	10/13	1
Euro OTC	Minoxidil	1807018	1
Euro OTC	Minoxidil	1102010-02	1
Euro OTC	Minoxidil	1612011-011	1
Euro OTC	Minoxidil	1112024-02	6
Euro OTC	Minoxidil	1612011	1
Euro OTC	Minoxidil	12/16	2
Euro OTC	Minoxidil	1112024-03	6
Euro OTC	Minoxidil	1612011-01	54
Euro OTC	Minoxidil	161201-01	1
Euro OTC	Minoxidil	1607001-01	1
Euro OTC	Minoxidil	140502801	4
Euro OTC	Minoxidil	16070007-01	1
Euro OTC	Minoxidil	1504013-01	16
Euro OTC	Minoxidil	700158	1
Euro OTC	Minoxidil	1707018-01	1
Euro OTC	Minoxidil	1405028-01	34
Euro OTC	Minoxidil	1509035-01	32
Euro OTC	Minoxidil	-B.:1607007-01	1
Euro OTC	Minoxidil	19/16	1
Euro OTC	Minoxidil	160700701	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Minoxidil	180708-01	1
Euro OTC	Minoxidil	1504013-02	26
Euro OTC	Minoxidil	1705008-01	41
Euro OTC	Minoxidil	1112024-01	2
Euro OTC	Minoxidil	1509035-01	3
Euro OTC	Minoxidil	1803004-01	23
Euro OTC	Minoxidil	6614849	1
Euro OTC	Minoxidil	1304029-01	18
Euro OTC	Minoxidil	1009042-01	3
Euro OTC	Minoxidil	1710043_01,180509M	1
Euro OTC	Minoxidil	14050281	1
Euro OTC	Minoxidil	2002312844	1
Euro OTC	Minoxidil	13001	1
Euro OTC	Minoxidil	17050008	1
Euro OTC	Minoxidil	171004301	5
Euro OTC	Minoxidil	16120-01	1
Euro OTC	Minoxidil	17100-43-01	1
Euro OTC	Minoxidil	13110838-01	1
Euro OTC	Minoxidil	1405028-01Mg:10g	1
Euro OTC	Minoxidil	130402901	2
Euro OTC	Minoxidil	141100601	5
Euro OTC	Minoxidil	12100056-01	1
Euro OTC	Minoxidil	1304029-C2	1
Euro OTC	Minoxidil	4/13	3
Euro OTC	Minoxidil	1405028-1	1
Euro OTC	Minoxidil	121005601	1
Euro OTC	Minoxidil	130402902	2
Euro OTC	Minoxidil	1210056-02	33
Euro OTC	Minoxidil	1602013-01	42
Euro OTC	Minoxidil	1411006601	1
Euro OTC	Minoxidil	1108039-03	1
Euro OTC	Minoxidil	1710043-01	45
Euro OTC	Mometasone furoate, m...	1203024-02	1
Euro OTC	Mometasone furoate, m...	1711027-01	7
Euro OTC	Mometasone furoate, m...	16j24-B11-332111	1
Euro OTC	Mometasone furoate, m...	161000901	1
Euro OTC	Mometasone furoate, m...	1711027	1
Euro OTC	Mometasone furoate, m...	1711027-02	1
Euro OTC	Mometasone furoate, m...	1810008-01	4
Euro OTC	Mometasone furoate, m...	17D19B02347349	1
Euro OTC	Mometasone furoate, m...	1610009-01	14
Euro OTC	Mometasone furoate, m...	160602101	1
Euro OTC	Mometasone furoate, m...	1606021-01	5
Euro OTC	Mometasone furoate, m...	1504012-02	1
Euro OTC	Mometasone furoate, m...	1510036-01	8
Euro OTC	Mometasone furoate, m...	1504012-01	6
Euro OTC	Mometasone furoate, m...	1408010-01	2
Euro OTC	Mometasone furoate, m...	5991	1
Euro OTC	Naphazoline hydrochlo...	1506015-02	1
Euro OTC	Naphazoline hydrochlo...	1307020-03	1
Euro OTC	Naphazoline hydrochlo...	616620001	1
Euro OTC	Naphazoline hydrochlo...	1407034-01	2
Euro OTC	Naphazoline hydrochlo...	1407034-02	4
Euro OTC	Naphazoline hydrochlo...	1607041-01	6
Euro OTC	Naphazoline hydrochlo...	1307020-03	1
Euro OTC	Naphazoline hydrochlo...	14070034-01	1
Euro OTC	Naphazoline hydrochlo...	1607041	1
Euro OTC	Naphazoline hydrochlo...	1506015-01	2
Euro OTC	Naphazoline hydrochlo...	150601502	1
Euro OTC	Naphazoline hydrochlo...	113031-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Naphazoline hydrochlo...	1206037-02	4
Euro OTC	Naphazoline hydrochlo...	1103031-02	1
Euro OTC	Naphazoline hydrochlo...	1506015-02	6
Euro OTC	Naphazoline hydrochlo...	1407034-03	1
Euro OTC	Naphazoline hydrochlo...	1206037-01	1
Euro OTC	Neomycin sulfate	1502004	1
Euro OTC	Neomycin sulfate	1102043-01	3
Euro OTC	Neomycin sulfate	1502004-03	20
Euro OTC	Neomycin sulfate	1102043-02	13
Euro OTC	Neomycin sulfate	1102043	2
Euro OTC	Neomycin sulfate	G73890	1
Euro OTC	Neomycin sulfate	1303033-03	3
Euro OTC	Neomycin sulfate	16019046-01	1
Euro OTC	Neomycin sulfate	N05780	2
Euro OTC	Neomycin sulfate	1303033-02	3
Euro OTC	Neomycin sulfate	1502004-02	5
Euro OTC	Neomycin sulfate	160104601	2
Euro OTC	Neomycin sulfate	1303033-01	12
Euro OTC	Neomycin sulfate	1601046-02	29
Euro OTC	Neomycin sulfate	151116N	1
Euro OTC	Neomycin sulfate	1601046-01	30
Euro OTC	Neomycin sulfate	L34419	3
Euro OTC	Neomycin sulfate	1601046-1	1
Euro OTC	Neomycin sulfate	1502004-01	17
Euro OTC	Neomycin sulfate	1802005-01	23
Euro OTC	Neomycin sulfate	1303001-03	1
Euro OTC	Neomycin sulfate	1303033-03	16
Euro OTC	Neomycin sulfate	160146-02	1
Euro OTC	Nettle seed	12131507	2
Euro OTC	Nicotinamide	1211031-02	6
Euro OTC	Nicotinamide	11208021	1
Euro OTC	Nicotinamide	1610019-01	8
Euro OTC	Nicotinamide	2610019-01	1
Euro OTC	Nicotinamide	1601045-01	4
Euro OTC	Nicotinamide	1108021-01	1
Euro OTC	Nicotinamide	1707001-ä01	1
Euro OTC	Nicotinamide	131002501	1
Euro OTC	Nicotinamide	1310025-01	3
Euro OTC	Nicotinamide	1801039-01	5
Euro OTC	Nicotinamide	1707001-01	9
Euro OTC	Nicotinamide	11401052	1
Euro OTC	Nicotinamide	1408019-19	1
Euro OTC	Nicotinamide	1503024-02	1
Euro OTC	Norethisterone acetate	1306033-02	2
Euro OTC	Nystatin	C129854	1
Euro OTC	Nystatin	1511020-02	1
Euro OTC	Nystatin	1307053-01	6
Euro OTC	Nystatin	1712010-01	23
Euro OTC	Nystatin	1107046-01	6
Euro OTC	Nystatin	161200103	2
Euro OTC	Nystatin	4014991	1
Euro OTC	Nystatin	16120001-02	1
Euro OTC	Nystatin	130705301	1
Euro OTC	Nystatin	121024-01	1
Euro OTC	Nystatin	171201001	1
Euro OTC	Nystatin	1511020-01	1
Euro OTC	Nystatin	1712010	1
Euro OTC	Nystatin	1712010-02	5
Euro OTC	Nystatin	171201002	3
Euro OTC	Nystatin	C129531	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Nystatin	1212024-01	15
Euro OTC	Nystatin	1207009-02	2
Euro OTC	Nystatin	19500711	1
Euro OTC	Nystatin	127009-02	1
Euro OTC	Nystatin	2811032-01	1
Euro OTC	Nystatin	1107036-01	6
Euro OTC	Nystatin	1307053-02	3
Euro OTC	Nystatin	C129440	1
Euro OTC	Nystatin	1212024-02	15
Euro OTC	Oak bark, cut		1
Euro OTC	Ornithine hydrochloride	516Q-03475	1
Euro OTC	Ornithine hydrochloride	4712A-03475	1
Euro OTC	Ornithine hydrochloride	1014I-03475	1
Euro OTC	Ornithine hydrochloride	1704006-01	1
Euro OTC	Ornithine hydrochloride	315E-03475	1
Euro OTC	Oxytetracycline hydro...	1801033-02	5
Euro OTC	Oxytetracycline hydro...	120502802	2
Euro OTC	Oxytetracycline hydro...	YT150502056	1
Euro OTC	Oxytetracycline hydro...	170034-01	1
Euro OTC	Oxytetracycline hydro...	1701034-015415	1
Euro OTC	Oxytetracycline hydro...	170103401	2
Euro OTC	Oxytetracycline hydro...	1407012-09	1
Euro OTC	Oxytetracycline hydro...	1407012-02	14
Euro OTC	Oxytetracycline hydro...	28011505	1
Euro OTC	Oxytetracycline hydro...	170103402	2
Euro OTC	Oxytetracycline hydro...	181016183	1
Euro OTC	Oxytetracycline hydro...	1801033-01	25
Euro OTC	Oxytetracycline hydro...	180103301	3
Euro OTC	Oxytetracycline hydro...	1306022-02	17
Euro OTC	Oxytetracycline hydro...	130602201	1
Euro OTC	Oxytetracycline hydro...	1205028	1
Euro OTC	Oxytetracycline hydro...	1306022	1
Euro OTC	Oxytetracycline hydro...	YT160702059	2
Euro OTC	Oxytetracycline hydro...	170612181	1
Euro OTC	Oxytetracycline hydro...	1701034-01	39
Euro OTC	Oxytetracycline hydro...	1707010TH, 1701034-01	1
Euro OTC	Oxytetracycline hydro...	15060109	1
Euro OTC	Oxytetracycline hydro...	15060377-01	1
Euro OTC	Oxytetracycline hydro...	150603701	1
Euro OTC	Oxytetracycline hydro...	1506037-01	29
Euro OTC	Oxytetracycline hydro...	180103302	1
Euro OTC	Oxytetracycline hydro...	1306022-03	5
Euro OTC	Oxytetracycline hydro...	1407012-03	6
Euro OTC	Oxytetracycline hydro...	1801033	1
Euro OTC	Oxytetracycline hydro...	15051507	1
Euro OTC	Oxytetracycline hydro...	1506037-03	18
Euro OTC	Oxytetracycline hydro...	1602025-01	19
Euro OTC	Oxytetracycline hydro...	160202501	2
Euro OTC	Oxytetracycline hydro...	1602025-02	51
Euro OTC	Oxytetracycline hydro...	1703130H/1602025-01	1
Euro OTC	Oxytetracycline hydro...	1602025	3
Euro OTC	Oxytetracycline hydro...	1602025-1	1
Euro OTC	Oxytetracycline hydro...	1407011202	1
Euro OTC	Oxytetracycline hydro...	375-18	1
Euro OTC	Oxytetracycline hydro...	1506037-01	3
Euro OTC	Oxytetracycline hydro...	1205028-01	10
Euro OTC	Oxytetracycline hydro...	1306022-01	25
Euro OTC	Oxytetracycline hydro...	1011045-01	3
Euro OTC	Oxytetracycline hydro...	1111061-01	7
Euro OTC	Oxytetracycline hydro...	12050-28	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Oxytetracycline hydro...	513E-02865	1
Euro OTC	Oxytetracycline hydro...	1111061-61	1
Euro OTC	Oxytetracycline hydro...	21111815	1
Euro OTC	Paracetamol	1111024-01	1
Euro OTC	Paracetamol	1111024-02	1
Euro OTC	Paracetamol	17028409	3
Euro OTC	Pheniramine maleate	1411044-01	1
Euro OTC	Pheniramine maleate	1411044-02	3
Euro OTC	Pheniramine maleate	1701021-01	6
Euro OTC	Pheniramine maleate	1209032-01	2
Euro OTC	Phenylbutazon	19091806	1
Euro OTC	Phenylbutazon	170612PB/1604033-01	1
Euro OTC	Pilocarpine hydrochlo...	Euroo	1
Euro OTC	Pilocarpine hydrochlo...	1603047-02	21
Euro OTC	Pilocarpine hydrochlo...	1806012-01	3
Euro OTC	Pilocarpine hydrochlo...	1603047	1
Euro OTC	Pilocarpine hydrochlo...	160304702	1
Euro OTC	Pilocarpine hydrochlo...	1603047-01	4
Euro OTC	Pilocarpine hydrochlo...	1305044-03	3
Euro OTC	Pilocarpine hydrochlo...	13050-14-01	1
Euro OTC	Pilocarpine hydrochlo...	130504401	3
Euro OTC	Pilocarpine hydrochlo...	1202031-01	1
Euro OTC	Pilocarpine hydrochlo...	1305044-01	4
Euro OTC	Pilocarpine hydrochlo...	1202031-02	1
Euro OTC	Polymyxin B sulfate	140601202	3
Euro OTC	Polymyxin B sulfate	28041730	1
Euro OTC	Polymyxin B sulfate	1604016-02	3
Euro OTC	Polymyxin B sulfate	1802006-01	6
Euro OTC	Polymyxin B sulfate	150500501	1
Euro OTC	Polymyxin B sulfate	1505005-01	4
Euro OTC	Polymyxin B sulfate	1304033-03	2
Euro OTC	Polymyxin B sulfate	23051502	1
Euro OTC	Polymyxin B sulfate	1406012-02	4
Euro OTC	Polymyxin B sulfate	1406012-01	1
Euro OTC	Polymyxin B sulfate	1406012-01	2
Euro OTC	Polymyxin B sulfate	180200601	1
Euro OTC	Polymyxin B sulfate	1304033-02	2
Euro OTC	Polymyxin B sulfate	11091419	1
Euro OTC	Polymyxin B sulfate	25101201	1
Euro OTC	Polymyxin B sulfate	1202039-01	3
Euro OTC	Polymyxin B sulfate	1701010-01	13
Euro OTC	Polymyxin B sulfate	1304033-01	1
Euro OTC	Potassium citrate	1406025-01	1
Euro OTC	Potassium citrate	1915A-01737	1
Euro OTC	Potassium citrate	1504024-01	2
Euro OTC	Potassium citrate	1311036-01	10
Euro OTC	Potassium citrate	22101303	1
Euro OTC	Potassium citrate	160604-01	1
Euro OTC	Potassium citrate	123456	7
Euro OTC	Potassium citrate	1603040-01	12
Euro OTC	Potassium citrate	1611017-01	5
Euro OTC	Potassium citrate	1705024-01	5
Euro OTC	Potassium citrate	1410014-01	2
Euro OTC	Potassium citrate	1306032-01	1
Euro OTC	Potassium citrate	1509044-01	5
Euro OTC	Potassium citrate	1606041-01	4
Euro OTC	Potassium citrate	160601-01	1
Euro OTC	Potassium citrate	160604101	2
Euro OTC	Potassium hydrocarbon...	1204038-01	2
Euro OTC	Potassium hydrocarbon...	1202013-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Potassium hydrocarbon...	2809033-01	1
Euro OTC	Potassium hydrocarbon...	1208010-01	3
Euro OTC	Potassium hydrocarbon...	4413I1728	1
Euro OTC	Potassium sorbate	1710012-01	1
Euro OTC	Potassium sorbate	1608013-01	5
Euro OTC	Potassium sorbate	216331021	8
Euro OTC	Potassium sorbate	180529-01	1
Euro OTC	Potassium sorbate	1111038-02	2
Euro OTC	Potassium sorbate	111103802	1
Euro OTC	Potassium sorbate	1312024-01	5
Euro OTC	Potassium sorbate	1508025-01	7
Euro OTC	Potassium sorbate	1111038-01	1
Euro OTC	Povidone iodine	8051808	1
Euro OTC	Povidone iodine	1110013-02	8
Euro OTC	Povidone iodine	1509009-01	1
Euro OTC	Povidone iodine	1704007-01	1
Euro OTC	Povidone iodine	1509009-02	2
Euro OTC	Povidone iodine	1606026-02	4
Euro OTC	Povidone iodine	1407032-01	6
Euro OTC	Povidone iodine	111001302	1
Euro OTC	Povidone iodine	1401029-01	1
Euro OTC	Povidone iodine	1805007-01	1
Euro OTC	Povidone iodine	1603023-01	2
Euro OTC	Povidone iodine	1606026	1
Euro OTC	Povidone iodine	G4208BPZ0	2
Euro OTC	Povidone iodine	130901001	1
Euro OTC	Povidone iodine	110013-02	1
Euro OTC	Prasterone	T03-04/1412028-01	1
Euro OTC	Prasterone	18050118-01	1
Euro OTC	Prasterone	180901001	1
Euro OTC	Prasterone	2081302	1
Euro OTC	Prasterone	1306021-01	3
Euro OTC	Prasterone	9215-4/09	1
Euro OTC	Prasterone	5101849	1
Euro OTC	Prasterone	23021711	1
Euro OTC	Prasterone	1602018-02	2
Euro OTC	Prasterone	1311002-01	2
Euro OTC	Prasterone	ZDHAC18002	1
Euro OTC	Prasterone	170202501	1
Euro OTC	Prasterone	1409004-01	2
Euro OTC	Prasterone	1401040-01	1
Euro OTC	Prasterone	XPC0263	1
Euro OTC	Prasterone	30041509	1
Euro OTC	Prasterone	26031506	1
Euro OTC	Prasterone	22021306	1
Euro OTC	Prasterone	DHEA150801-1	1
Euro OTC	Prasterone	1602018-01	3
Euro OTC	Prasterone	1702025-01	5
Euro OTC	Prasterone	1702025	1
Euro OTC	Prasterone	1704012-01	1
Euro OTC	Prasterone	1506011-01	3
Euro OTC	Prasterone	4091511	1
Euro OTC	Prasterone	1412028	1
Euro OTC	Prasterone	1805018-01	2
Euro OTC	Prasterone	1809010-01	1
Euro OTC	Prasterone	1412028-01	5
Euro OTC	Prasterone	1203010-01	1
Euro OTC	Prasterone	1409004-01	2
Euro OTC	Prednicarbate, micron...	1809015-01	11
Euro OTC	Prednicarbate, micron...	180901501	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Prednicarbate, micron...	170702201	1
Euro OTC	Prednicarbate, micron...	180101301	1
Euro OTC	Prednicarbate, micron...	1801013	1
Euro OTC	Prednicarbate, micron...	1805001-01	2
Euro OTC	Prednicarbate, micron...	1805005-01	9
Euro OTC	Prednicarbate, micron...	180500501	1
Euro OTC	Prednicarbate, micron...	1707022-01	27
Euro OTC	Prednicarbate, micron...	1501033-02	4
Euro OTC	Prednicarbate, micron...	1707022-01, 180425PC	1
Euro OTC	Prednicarbate, micron...	1801013-01	30
Euro OTC	Prednicarbate, micron...	1302033-02	2
Euro OTC	Prednicarbate, micron...	1302033-01	1
Euro OTC	Prednicarbate, micron...	1501033-01	2
Euro OTC	Prednicarbate, micron...	1507011-02	4
Euro OTC	Prednicarbate, micron...	1601058-01	4
Euro OTC	Prednicarbate, micron...	1605001	1
Euro OTC	Prednicarbate, micron...	1809015001	1
Euro OTC	Prednicarbate, micron...	1903023-01	1
Euro OTC	Prednicarbate, micron...	160902	1
Euro OTC	Prednicarbate, micron...	13224402	1
Euro OTC	Prednicarbate, micron...	1707022201	1
Euro OTC	Prednicarbate, micron...	1605001-01	7
Euro OTC	Prednicarbate, micron...	1609022-01	23
Euro OTC	Prednicarbate, micron...	1609011-01	1
Euro OTC	Prednisolone acetate	1306040-01	23
Euro OTC	Prednisolone acetate	1701041-02	19
Euro OTC	Prednisolone acetate	170104102	2
Euro OTC	Prednisolone acetate	1805009-01	23
Euro OTC	Prednisolone acetate	18004403	2
Euro OTC	Prednisolone acetate	1306040-1	1
Euro OTC	Prednisolone acetate	1701041-03	4
Euro OTC	Prednisolone acetate	1805009-02	10
Euro OTC	Prednisolone acetate	1701041-01, 180410PA	1
Euro OTC	Prednisolone acetate	2(1903031-01)	1
Euro OTC	Prednisolone acetate	1805009	2
Euro OTC	Prednisolone acetate	1209028	1
Euro OTC	Prednisolone acetate	1903031-01	3
Euro OTC	Prednisolone acetate	514U2106	1
Euro OTC	Prednisolone acetate	1209028-02	18
Euro OTC	Prednisolone acetate	200540	2
Euro OTC	Prednisolone acetate	18067320	1
Euro OTC	Prednisolone acetate	K07D20161003M	1
Euro OTC	Prednisolone acetate	2413Q02106	1
Euro OTC	Prednisolone acetate	190303101	1
Euro OTC	Prednisolone acetate	170413PA/1607010-01	1
Euro OTC	Prednisolone acetate	1501036-01	19
Euro OTC	Prednisolone acetate	1607010-10	1
Euro OTC	Prednisolone acetate	1406014-01	20
Euro OTC	Prednisolone acetate	1406012-02	1
Euro OTC	Prednisolone acetate	289295	1
Euro OTC	Prednisolone acetate	150103601	3
Euro OTC	Prednisolone acetate	2714U-02106	1
Euro OTC	Prednisolone acetate	1460614-01	1
Euro OTC	Prednisolone acetate	1203028-03	1
Euro OTC	Prednisolone acetate	1501036-02	31
Euro OTC	Prednisolone acetate	1501036	1
Euro OTC	Prednisolone acetate	K07A20140304M	1
Euro OTC	Prednisolone acetate	1701041-01	29
Euro OTC	Prednisolone acetate	5012I-02106	1
Euro OTC	Prednisolone acetate	1501036.02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Prednisolone acetate	1209028-03	2
Euro OTC	Prednisolone acetate	130604002	1
Euro OTC	Prednisolone acetate	150103602	1
Euro OTC	Prednisolone acetate	161005PA	1
Euro OTC	Prednisolone acetate	1607010-01	87
Euro OTC	Prednisolone acetate	2215I-02106	1
Euro OTC	Prednisolone acetate	1209028-01	5
Euro OTC	Prednisolone acetate	2507212106	1
Euro OTC	Prednisolone acetate	1108043-02	10
Euro OTC	Prednisolone acetate	711212106	1
Euro OTC	Prednisolone acetate	151202001	2
Euro OTC	Prednisolone acetate	1512020-01	35
Euro OTC	Prednisolone acetate	1108043-01	2
Euro OTC	Prednisolone, microni...	16E09-B07-329686	3
Euro OTC	Prednisolone, microni...	14/10-B03-302263	1
Euro OTC	Prednisolone, microni...	170202201	1
Euro OTC	Prednisolone, microni...	1702020-01	30
Euro OTC	Prednisolone, microni...	1702020-02	15
Euro OTC	Prednisolone, microni...	1501047-03	2
Euro OTC	Prednisolone, microni...	1702020-03	2
Euro OTC	Prednisolone, microni...	16134308	9
Euro OTC	Prednisolone, microni...	170202002	3
Euro OTC	Prednisolone, microni...	150104743	1
Euro OTC	Prednisolone, microni...	121221-03	1
Euro OTC	Prednisolone, microni...	1108055-01	1
Euro OTC	Prednisolone, microni...	1108055-03	1
Euro OTC	Prednisolone, microni...	110805504	3
Euro OTC	Prednisolone, microni...	1501047-01	28
Euro OTC	Prednisolone, microni...	150104701	5
Euro OTC	Prednisolone, microni...	1212021-03	36
Euro OTC	Prednisolone, microni...	1212021	1
Euro OTC	Prednisolone, microni...	1607011-01	30
Euro OTC	Prednisolone, microni...	15010477-01	1
Euro OTC	Prednisolone, microni...	1501047	1
Euro OTC	Prednisolone, microni...	501047-01	1
Euro OTC	Prednisolone, microni...	1501047-02	17
Euro OTC	Prednisolone, microni...	150104747-01	1
Euro OTC	Prednisolone, microni...	1212021-01	3
Euro OTC	Prednisolone, microni...	12021-03	1
Euro OTC	Prednisolone, microni...	1212021-02	4
Euro OTC	Prednisolone, microni...	121202101	1
Euro OTC	Prednisolone, microni...	121202102	1
Euro OTC	Prednisolone, microni...	121021-01	1
Euro OTC	Prednisolone, microni...	1212021-01	11
Euro OTC	Prednisolone, microni...	1108055-04	2
Euro OTC	Prednisolone, microni...	30042014B	1
Euro OTC	Prednisolone, microni...	1314Q2105	1
Euro OTC	Prednisolone, microni...	118055-04	1
Euro OTC	Prednisolone, microni...	NEPL110103	2
Euro OTC	Prednisolone, microni...	1212021-02	24
Euro OTC	Prednisolone, microni...	1108055-02	2
Euro OTC	Prednisone	1702008-01	1
Euro OTC	Prednisone	150043-03	1
Euro OTC	Prednisone	1401042-02	1
Euro OTC	Prednisone	1604051-01	1
Euro OTC	Prednisone	1501043-03	4
Euro OTC	Prednisone	1501043-02	1
Euro OTC	Prednisone	1604501-01	1
Euro OTC	Pregnenolone	1312007-01	1
Euro OTC	Pregnenolone	1402025-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Pregnenolone	1606028-01	10
Euro OTC	Pregnenolone	1602001-01	1
Euro OTC	Pregnenolone	1612019-01	51
Euro OTC	Procaine hydrochloride	1601014-02	1
Euro OTC	Procaine hydrochloride	1601014-02	10
Euro OTC	Procaine hydrochloride	200644	1
Euro OTC	Procaine hydrochloride	1104012-03	4
Euro OTC	Procaine hydrochloride	1003001-01	1
Euro OTC	Procaine hydrochloride	1104012-04	3
Euro OTC	Procaine hydrochloride	1104012-05	5
Euro OTC	Procaine hydrochloride	5127	1
Euro OTC	Procaine hydrochloride	1405027-02	1
Euro OTC	Procaine hydrochloride	1601014-01	8
Euro OTC	Procaine hydrochloride	1104012-01	1
Euro OTC	Procaine hydrochloride	1405027-01	4
Euro OTC	Progesterone, microni...	160501303	2
Euro OTC	Progesterone, microni...	170502801	1
Euro OTC	Progesterone, microni...	1406022-01	2
Euro OTC	Progesterone, microni...	1801036-01	23
Euro OTC	Progesterone, microni...	1705028	1
Euro OTC	Progesterone, microni...	1801036-02	10
Euro OTC	Progesterone, microni...	1806020-01	54
Euro OTC	Progesterone, microni...	1809017-01	55
Euro OTC	Progesterone, microni...	180901701	20
Euro OTC	Progesterone, microni...	1811033-01	64
Euro OTC	Progesterone, microni...	1901022-01	31
Euro OTC	Progesterone, microni...	1801036.01	1
Euro OTC	Progesterone, microni...	1706009-01	9
Euro OTC	Progesterone, microni...	1605013-03	4
Euro OTC	Progesterone, microni...	1705028-01	62
Euro OTC	Progesterone, microni...	1704012-01	8
Euro OTC	Progesterone, microni...	140910-01	1
Euro OTC	Progesterone, microni...	5104	1
Euro OTC	Progesterone, microni...	1111053-01	6
Euro OTC	Progesterone, microni...	1206024-01	1
Euro OTC	Progesterone, microni...	1211004-01	5
Euro OTC	Progesterone, microni...	co1-111218(M)	1
Euro OTC	Progesterone, microni...	1304040-01	5
Euro OTC	Progesterone, microni...	C01-130301M	1
Euro OTC	Progesterone, microni...	5102	1
Euro OTC	Progesterone, microni...	140602202	1
Euro OTC	Progesterone, microni...	1304040-02	3
Euro OTC	Progesterone, microni...	130404001	1
Euro OTC	Progesterone, microni...	c01-130301	1
Euro OTC	Progesterone, microni...	c01-111218	1
Euro OTC	Progesterone, microni...	C01-130301 (M)	2
Euro OTC	Progesterone, microni...	1308015-01	4
Euro OTC	Progesterone, microni...	1308015-01	1
Euro OTC	Progesterone, microni...	7406022-01	1
Euro OTC	Progesterone, microni...	1403019-01	2
Euro OTC	Progesterone, microni...	1409029-1	1
Euro OTC	Progesterone, microni...	1709012-01	99
Euro OTC	Progesterone, microni...	1406910-01	1
Euro OTC	Progesterone, microni...	1411043-01	4
Euro OTC	Progesterone, microni...	1406022-02	2
Euro OTC	Progesterone, microni...	170103802	9
Euro OTC	Progesterone, microni...	1701038-02	33
Euro OTC	Progesterone, microni...	170138-01	20
Euro OTC	Progesterone, microni...	1701038-01	30
Euro OTC	Progesterone, microni...	1609001-01	44

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Progesterone, microni...	160100901	1
Euro OTC	Progesterone, microni...	1605013-01	22
Euro OTC	Progesterone, microni...	HTT160102	1
Euro OTC	Progesterone, microni...	1601009-02	4
Euro OTC	Progesterone, microni...	1503033-01	1
Euro OTC	Progesterone, microni...	1601009-01	8
Euro OTC	Progesterone, microni...	1509008_01	1
Euro OTC	Progesterone, microni...	1509008-02	2
Euro OTC	Progesterone, microni...	1509008-01	3
Euro OTC	Progesterone, microni...	1506033-01	8
Euro OTC	Progesterone, microni...	1503012-01	6
Euro OTC	Progesterone, microni...	141104301	1
Euro OTC	Progesterone, microni...	1502030-01	2
Euro OTC	Progesterone, microni...	1409019-01	1
Euro OTC	Progesterone, microni...	1409029-01	3
Euro OTC	Progesterone, microni...	1605013-02	44
Euro OTC	Propranolol hydrochlo...	13379806	1
Euro OTC	Propranolol hydrochlo...	13I13-N04	1
Euro OTC	Propranolol hydrochlo...	16138304	12
Euro OTC	Propranolol hydrochlo...	15K11B01-320607	1
Euro OTC	Propyl-4-hydroxybenzo...	2509025	1
Euro OTC	Propyphenazone	B:09369110	1
Euro OTC	Pyridoxal-5-phosphate	1808030-02	1
Euro OTC	Pyridoxal-5-phosphate	1711030-01	2
Euro OTC	Pyridoxal-5-phosphate	1701015-01	3
Euro OTC	Pyridoxine hydrochlor...	1601039-01	1
Euro OTC	Pyridoxine hydrochlor...	160803501	1
Euro OTC	Pyridoxine hydrochlor...	1710014-01	1
Euro OTC	Pyridoxine hydrochlor...	1110001-01	1
Euro OTC	Pyridoxine hydrochlor...	17001401	1
Euro OTC	Pyridoxine hydrochlor...	1206022-02	1
Euro OTC	Pyridoxine hydrochlor...	1604063-01	1
Euro OTC	Riboflavin	1701050-02	1
Euro OTC	Riboflavin	L1210048-01	1
Euro OTC	Riboflavin	L1305039-01	1
Euro OTC	Riboflavin	1707017-01	1
Euro OTC	Salicylic acid	1303020-01	28
Euro OTC	Salicylic acid	1110020	1
Euro OTC	Salicylic acid	115465478596	1
Euro OTC	Salicylic acid	140602702	2
Euro OTC	Salicylic acid	1707009-01	17
Euro OTC	Salicylic acid	160502301	2
Euro OTC	Salicylic acid	160502331	1
Euro OTC	Salicylic acid	1605023-02	14
Euro OTC	Salicylic acid	1605023/02	1
Euro OTC	Salicylic acid	1511023B	1
Euro OTC	Salicylic acid	1110020-01	4
Euro OTC	Salicylic acid	16E02-B09-325356	1
Euro OTC	Salicylic acid	1406027-02	13
Euro OTC	Salicylic acid	1605023-01	19
Euro OTC	Salicylic acid	1509007/02	1
Euro OTC	Salicylic acid	170700901	3
Euro OTC	Salicylic acid	1406027-02Mg1kg	1
Euro OTC	Salicylic acid	1901015-01	1
Euro OTC	Salicylic acid	1707009-02	8
Euro OTC	Salicylic acid	1509007-01	28
Euro OTC	Salicylic acid	1408027	1
Euro OTC	Salicylic acid	1801042-01	16
Euro OTC	Salicylic acid	102455	1
Euro OTC	Salicylic acid	170700902	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Salicylic acid	1406027-01	11
Euro OTC	Salicylic acid	1303020-01Mg:1kg	1
Euro OTC	Salicylic acid	RAS1703500	3
Euro OTC	Salicylic acid	130302001	1
Euro OTC	Salicylic acid	1303020-02	4
Euro OTC	Serine	1606020-01	3
Euro OTC	Sha Ren	H0939022	1
Euro OTC	Sildenafil citrate	1311010-03	6
Euro OTC	Sildenafil citrate	13110-01	1
Euro OTC	Sildenafil citrate	1807037-01	1
Euro OTC	Sildenafil citrate	19011907	1
Euro OTC	Sildenafil citrate	14051402	1
Euro OTC	Sildenafil citrate	1311010-01	6
Euro OTC	Sildenafil citrate	1807037-01x2	1
Euro OTC	Sildenafil citrate	17061405	1
Euro OTC	Sildenafil citrate	17081504	1
Euro OTC	Sildenafil citrate	131101003	1
Euro OTC	Sildenafil citrate	311010-06	1
Euro OTC	Sildenafil citrate	1311010-6	1
Euro OTC	Sildenafil citrate	15081830	1
Euro OTC	Sildenafil citrate	103031328	1
Euro OTC	Sildenafil citrate	1311010	1
Euro OTC	Sildenafil citrate	27021802	1
Euro OTC	Sildenafil citrate	131101004	1
Euro OTC	Sildenafil citrate	1311010-06	8
Euro OTC	Sildenafil citrate	1311010-05	5
Euro OTC	Sildenafil citrate	10061706	1
Euro OTC	Sodium acetate trihyd...	1005020V-01	1
Euro OTC	Sodium alginate	1509036-01	1
Euro OTC	Sodium ascorbate	28081304	1
Euro OTC	Sodium ascorbate	11121301	1
Euro OTC	Sodium ascorbate	1210026-01	1
Euro OTC	Sodium ascorbate	7011405	1
Euro OTC	Sodium ascorbate	L1308005-01	1
Euro OTC	Sodium ascorbate	16091301	1
Euro OTC	Sodium ascorbate	4021401	1
Euro OTC	Sodium ascorbate	20091303	1
Euro OTC	Sodium ascorbate	13031402	1
Euro OTC	Sodium ascorbate	13051406	1
Euro OTC	Sodium ascorbate	12061407	1
Euro OTC	Sodium ascorbate	11081401	1
Euro OTC	Sodium ascorbate	170105401	1
Euro OTC	Sodium benzoate	1809005-01	1
Euro OTC	Sodium benzoate	1803034-01	1
Euro OTC	Sodium benzoate	1112006-01	2
Euro OTC	Sodium benzoate	171101901	1
Euro OTC	Sodium benzoate	170702401	1
Euro OTC	Sodium benzoate	1209018-18	1
Euro OTC	Sodium benzoate	180900501	1
Euro OTC	Sodium benzoate	1007019-01	1
Euro OTC	Sodium benzoate	1701045	1
Euro OTC	Sodium benzoate	130402401	2
Euro OTC	Sodium benzoate	1602020-01	3
Euro OTC	Sodium benzoate	15060119-01	1
Euro OTC	Sodium benzoate	1304024-01	2
Euro OTC	Sodium benzoate	1409040-01	1
Euro OTC	Sodium benzoate	151004401	3
Euro OTC	Sodium benzoate	1209018-01	5
Euro OTC	Sodium benzoate	1510044-01	11
Euro OTC	Sodium benzoate	1404011-01	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Sodium benzoate	2604040V-05	1
Euro OTC	Sodium benzoate	1307013-01	1
Euro OTC	Sodium benzoate	1309027-01	2
Euro OTC	Sodium benzoate	1607005-01	4
Euro OTC	Sodium benzoate	1307013-01	1
Euro OTC	Sodium citrate	1404033-01	6
Euro OTC	Sodium citrate	17061407	1
Euro OTC	Sodium citrate	2114M1900	1
Euro OTC	Sodium citrate	1073240	1
Euro OTC	Sodium citrate	1107011-01	2
Euro OTC	Sodium citrate	1311035-01	19
Euro OTC	Sodium citrate	30111301	1
Euro OTC	Sodium citrate	1305001-01	10
Euro OTC	Sodium citrate	2903001-01	1
Euro OTC	Sodium citrate	1211030-01	11
Euro OTC	Sodium citrate	1413M-01900	1
Euro OTC	Sodium citrate	20200810-6	1
Euro OTC	Sodium citrate	13I24-B13-287921	1
Euro OTC	Sodium citrate	1404033-01	2
Euro OTC	Sodium citrate	140403301	5
Euro OTC	Sodium citrate		7
Euro OTC	Sodium citrate	156422	1
Euro OTC	Sodium citrate	2546426	1
Euro OTC	Sodium citrate	3454	1
Euro OTC	Sodium citrate	1645265	3
Euro OTC	Sodium citrate	358528	1
Euro OTC	Sodium citrate	815757/6745	1
Euro OTC	Sodium citrate	25425	2
Euro OTC	Sodium citrate	1311025-01	1
Euro OTC	Sodium citrate	1521017-01	1
Euro OTC	Sodium citrate	15111017-01	1
Euro OTC	Sodium citrate	1607040-01	21
Euro OTC	Sodium citrate	160740-01	1
Euro OTC	Sodium citrate	345678	1
Euro OTC	Sodium citrate	16070-01	1
Euro OTC	Sodium citrate	3032805	2
Euro OTC	Sodium citrate	1610016-01	33
Euro OTC	Sodium citrate	1901035-01	1
Euro OTC	Sodium citrate	1807035-01ciclo	1
Euro OTC	Sodium citrate	1807035-02	1
Euro OTC	Sodium citrate	48(08598497)	1
Euro OTC	Sodium citrate	3048768	2
Euro OTC	Sodium citrate	1807035-01	12
Euro OTC	Sodium citrate	1311035-01	4
Euro OTC	Sodium citrate	1710006-02	4
Euro OTC	Sodium citrate	5116E-01900	1
Euro OTC	Sodium citrate	4716E-01900	1
Euro OTC	Sodium citrate	161001601	1
Euro OTC	Sodium citrate	1610016-01	1
Euro OTC	Sodium citrate	1701029-01	30
Euro OTC	Sodium citrate	1610016601	1
Euro OTC	Sodium citrate	1701029-02	1
Euro OTC	Sodium citrate	151101701	2
Euro OTC	Sodium citrate	12345	2
Euro OTC	Sodium citrate	22101305	1
Euro OTC	Sodium citrate	1511017-01	34
Euro OTC	Sodium citrate	1412014-01	1
Euro OTC	Sodium citrate	23062015C	1
Euro OTC	Sodium citrate	1305001-01	1
Euro OTC	Sodium citrate	1205026-02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Sodium citrate	1311035-01	2
Euro OTC	Sodium citrate	1506041-01	16
Euro OTC	Sodium citrate	14201401	1
Euro OTC	Sodium citrate	1305001-01	1
Euro OTC	Sodium citrate	1506041-01	1
Euro OTC	Sodium citrate	1506041-01	1
Euro OTC	Sodium cromoglicate	1403017-01	1
Euro OTC	Sodium cromoglicate	170302701	1
Euro OTC	Sodium dihydrogen pho...	1305028-01	2
Euro OTC	Sodium dihydrogen pho...	1209004-01	2
Euro OTC	Sodium hydrogen carbo...	1304017-02	2
Euro OTC	Sodium hydrogen carbo...	1012011-01	1
Euro OTC	Sodium hydrogen carbo...	1109029-01	8
Euro OTC	Sodium hydrogen carbo...	1104013-01	2
Euro OTC	Sodium hydrogen carbo...	120800801	1
Euro OTC	Sodium hydrogen carbo...	1304017-01	5
Euro OTC	Sodium hydrogen carbo...	1012011	1
Euro OTC	Sodium hydrogen carbo...	120303001	1
Euro OTC	Sodium hydrogen carbo...	1208008-01	9
Euro OTC	Sodium sulfate decahy...	1108035-02	2
Euro OTC	Sodium sulfate decahy...	1210005-02	2
Euro OTC	Sodium sulfate decahy...	1108035-01	2
Euro OTC	Sodium sulfate decahy...	4312I-01920	1
Euro OTC	Sodium sulfate decahy...	121000502	2
Euro OTC	Sodium tetraborate	1701046-01,170913NTB	1
Euro OTC	Sodium tetraborate	1711002-01	3
Euro OTC	Sodium tetraborate	1701046-01	1
Euro OTC	Sodium tetraborate	1607036-01	4
Euro OTC	Sodium tetraborate	1507024-01	3
Euro OTC	Sodium tetraborate	1309032-01	2
Euro OTC	Sodium tetraborate	130203401	1
Euro OTC	Sodium tetraborate	1302034-01	1
Euro OTC	Sodium tetraborate	160304301	1
Euro OTC	Sodium tetraborate	1408006-01	1
Euro OTC	Sodium tetraborate	1503043-01	1
Euro OTC	Sodium tetraborate	13111207	1
Euro OTC	Sorbic acid	o216E-01048	1
Euro OTC	Sorbitol	1610018-01	11
Euro OTC	Sorbitol	1610018-03	3
Euro OTC	Sorbitol	1604007-01	5
Euro OTC	Sorbitol	150652-01	1
Euro OTC	Sorbitol	1501048-01	3
Euro OTC	Sorbitol	1406014-01	1
Euro OTC	Sorbitol	700476	1
Euro OTC	Sorbitol	1301019-01	4
Euro OTC	Sorbitol	1311034-01	1
Euro OTC	Sorbitol	1404007-01	2
Euro OTC	Sorbitol	1506052-01	4
Euro OTC	Sorbitol	1610018-02	1
Euro OTC	Spironolactone	1212034-01	3
Euro OTC	Spironolactone	1809020-01	1
Euro OTC	Spironolactone	161200201	2
Euro OTC	Spironolactone	16020006-01	1
Euro OTC	Spironolactone	1711013-01	1
Euro OTC	Spironolactone	1000403502	1
Euro OTC	Spironolactone	1612002-01	5
Euro OTC	Spironolactone	3513E-03249	1
Euro OTC	Spironolactone	1101013-01	2
Euro OTC	Spironolactone	1504016-01	2
Euro OTC	Spironolactone	1402029-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Spironolactone	1004035	1
Euro OTC	Sucrose	12345	2
Euro OTC	Sulfacetamide	1312012-01	6
Euro OTC	Sulfacetamide	1312012-02	4
Euro OTC	Sulfacetamide	1509018-01	4
Euro OTC	Sulfacetamide	1509018-02	4
Euro OTC	Sulfacetamide	1702006-01	1
Euro OTC	Sulfacetamide	1702006-02	2
Euro OTC	Sulfacetamide	1808035-01	1
Euro OTC	Sulfacetamide	1112032-02	4
Euro OTC	Tannin	16074501	1
Euro OTC	Tannin	16074509	8
Euro OTC	Tartaric acid	1604046-02	1
Euro OTC	Tartaric acid	702265	1
Euro OTC	Taurine	1704001-01	3
Euro OTC	Taurine	XP16023196	1
Euro OTC	Testosterone	1106011-03	1
Euro OTC	Testosterone	1307001-01	6
Euro OTC	Testosterone	1207005-02	1
Euro OTC	Testosterone	1106011-01	1
Euro OTC	Testosterone	1106011-03	6
Euro OTC	Testosterone	1404002-01	1
Euro OTC	Testosterone	1607013-01	1
Euro OTC	Testosterone	1607013-02	6
Euro OTC	Testosterone	18147813	2
Euro OTC	Testosterone	1811027-01	4
Euro OTC	Testosterone	1808033-01	1
Euro OTC	Testosterone	1805021-01	4
Euro OTC	Testosterone	140400201	1
Euro OTC	Testosterone	1404002-1	1
Euro OTC	Testosterone	1805006-01	1
Euro OTC	Testosterone	1511005-01	3
Euro OTC	Testosterone	180506-01	1
Euro OTC	Testosterone	1604049-01	4
Euro OTC	Testosterone	1802015	1
Euro OTC	Testosterone	APZTSTG18003T	9
Euro OTC	Testosterone	1702017-01	39
Euro OTC	Testosterone	1404002-01	1
Euro OTC	Testosterone	3x1106011-03	1
Euro OTC	Testosterone propionate	1512023-01	2
Euro OTC	Testosterone propionate	970116/1	1
Euro OTC	Testosterone propionate	1109024-02	3
Euro OTC	Testosterone propionate	1109024-03	1
Euro OTC	Testosterone propionate	1409037-01	1
Euro OTC	Testosterone propionate	1103010-02	1
Euro OTC	Testosterone propionate	1504023-01	2
Euro OTC	Testosterone propionate	15444718	2
Euro OTC	Testosterone propionate	1511003-01	7
Euro OTC	Testosterone propionate	1409037-01	1
Euro OTC	Testosterone propionate	1103010-01	1
Euro OTC	Testosterone propionate	1811040-01	1
Euro OTC	Tetracaine base	1803002-01	1
Euro OTC	Tetracaine base	151022-01	1
Euro OTC	Tetracaine base	1811021-01	2
Euro OTC	Tetracaine base	18041809	1
Euro OTC	Tetracaine base	1708024-01	5
Euro OTC	Tetracaine base	12081501	1
Euro OTC	Tetracaine base	1312034-01	7
Euro OTC	Tetracaine base	3081701	1
Euro OTC	Tetracaine base	1701035-01	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Tetracaine base	TC0091115	1
Euro OTC	Tetracaine base	1405025-01	2
Euro OTC	Tetracaine base	160600-01	1
Euro OTC	Tetracaine base	160901	1
Euro OTC	Tetracaine base	1512022-01	3
Euro OTC	Tetracaine base	1609010-01	13
Euro OTC	Tetracaine base	10887133	1
Euro OTC	Tetracaine base	1606009-01	8
Euro OTC	Tetracaine hydrochlor...	1603026-01	1
Euro OTC	Tetracaine hydrochlor...	20160812-2	1
Euro OTC	Tetracaine hydrochlor...	1110010-02	8
Euro OTC	Tetracaine hydrochlor...	3071716	1
Euro OTC	Tetracaine hydrochlor...	160703101	2
Euro OTC	Tetracaine hydrochlor...	1712012-01	11
Euro OTC	Tetracaine hydrochlor...	1607032-01	1
Euro OTC	Tetracaine hydrochlor...	1405026-04	2
Euro OTC	Tetracaine hydrochlor...	130301201	2
Euro OTC	Tetracaine hydrochlor...	1607031	1
Euro OTC	Tetracaine hydrochlor...	1409026-02	15
Euro OTC	Tetracaine hydrochlor...	14090-26	1
Euro OTC	Tetracaine hydrochlor...	170323181	1
Euro OTC	Tetracaine hydrochlor...	1409026	1
Euro OTC	Tetracaine hydrochlor...	1612005-015445	1
Euro OTC	Tetracaine hydrochlor...	1411019-01	29
Euro OTC	Tetracaine hydrochlor...	1409026-01	13
Euro OTC	Tetracaine hydrochlor...	5016Y-02039	1
Euro OTC	Tetracaine hydrochlor...	11110010-02	1
Euro OTC	Tetracaine hydrochlor...	24102014A	1
Euro OTC	Tetracaine hydrochlor...	140502601	1
Euro OTC	Tetracaine hydrochlor...	1207001-01	11
Euro OTC	Tetracaine hydrochlor...	2906009V-03	1
Euro OTC	Tetracaine hydrochlor...	1303012-01	25
Euro OTC	Tetracaine hydrochlor...	1310006-01	23
Euro OTC	Tetracaine hydrochlor...	1405026-01	20
Euro OTC	Tetracaine hydrochlor...	1310006-02	1
Euro OTC	Tetracaine hydrochlor...	1607031-01	40
Euro OTC	Tetracaine hydrochlor...	10416THCL	1
Euro OTC	Tetracaine hydrochlor...	1411019	1
Euro OTC	Tetracaine hydrochlor...	15120019-01	1
Euro OTC	Tetracaine hydrochlor...	1506040-02	7
Euro OTC	Tetracaine hydrochlor...	1512019-01	27
Euro OTC	Tetracaine hydrochlor...	140902601	1
Euro OTC	Tetracaine hydrochlor...	141101901	1
Euro OTC	Tetracaine hydrochlor...	4811A-02039	1
Euro OTC	Tetracaine hydrochlor...	3113M-02039	1
Euro OTC	Tetracaine hydrochlor...	150604001	2
Euro OTC	Tetracaine hydrochlor...	1812006-01	5
Euro OTC	Tetracaine hydrochlor...	17081503	1
Euro OTC	Tetracaine hydrochlor...	181200601	1
Euro OTC	Tetracaine hydrochlor...	1506040-01	29
Euro OTC	Tetracaine hydrochlor...	13509155	1
Euro OTC	Tetracaine hydrochlor...	161200501	2
Euro OTC	Tetracaine hydrochlor...	1106035	1
Euro OTC	Tetracaine hydrochlor...	151-2019-01	1
Euro OTC	Tetracaine hydrochlor...	1612005-01	25
Euro OTC	Tetracaine hydrochlor...	1803001-1	1
Euro OTC	Tetracaine hydrochlor...	21071807	1
Euro OTC	Tetracaine hydrochlor...	1303012-01Mg:10g	1
Euro OTC	Tetracaine hydrochlor...	180300101	5
Euro OTC	Tetracaine hydrochlor...	1803001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Tetracaine hydrochlor...	14G14-B08-300202	1
Euro OTC	Tetracaine hydrochlor...	1207001-02	10
Euro OTC	Tetracaine hydrochlor...	18701701	1
Euro OTC	Tetracaine hydrochlor...	1405026-02	1
Euro OTC	Tetracaine hydrochlor...	13155710	4
Euro OTC	Tetracaine hydrochlor...	1803001-01	15
Euro OTC	Tetracaine hydrochlor...	1807017-01	3
Euro OTC	Tetracaine hydrochlor...	1512019-02	2
Euro OTC	Tetracycline hydrochl...	141200701	1
Euro OTC	Tetracycline hydrochl...	140410-62	1
Euro OTC	Tetracycline hydrochl...	1506012-01	31
Euro OTC	Tetracycline hydrochl...	160621TCH	1
Euro OTC	Tetracycline hydrochl...	150601201	2
Euro OTC	Tetracycline hydrochl...	1811031-01	5
Euro OTC	Tetracycline hydrochl...	1706020-02	1
Euro OTC	Tetracycline hydrochl...	22(1706020-01)	1
Euro OTC	Tetracycline hydrochl...	170602001	1
Euro OTC	Tetracycline hydrochl...	17060201-01	1
Euro OTC	Tetracycline hydrochl...	1706020-01	27
Euro OTC	Tetracycline hydrochl...	180306T/1412007-03	1
Euro OTC	Tetracycline hydrochl...	411682	1
Euro OTC	Tetracycline hydrochl...	160302	1
Euro OTC	Tetracycline hydrochl...	1506012-02Eur	1
Euro OTC	Tetracycline hydrochl...	27101708	1
Euro OTC	Tetracycline hydrochl...	1506012-02	3
Euro OTC	Tetracycline hydrochl...	1412007-03	8
Euro OTC	Tetracycline hydrochl...	160302601	1
Euro OTC	Tetracycline hydrochl...	21031702	1
Euro OTC	Tetracycline hydrochl...	1603026	2
Euro OTC	Tetracycline hydrochl...	1603026_01	1
Euro OTC	Tetracycline hydrochl...	13042015C	1
Euro OTC	Tetracycline hydrochl...	1603026-01	43
Euro OTC	Tetracycline hydrochl...	14120007-01	1
Euro OTC	Tetracycline hydrochl...	21111404	1
Euro OTC	Tetracycline hydrochl...	1412007-01	12
Euro OTC	Tetracycline hydrochl...	140203001	1
Euro OTC	Tetracycline hydrochl...	100518-62	1
Euro OTC	Tetracycline hydrochl...	1204039-01	6
Euro OTC	Tetracycline hydrochl...	1103003-01	1
Euro OTC	Tetracycline hydrochl...	12040390	1
Euro OTC	Tetracycline hydrochl...	28021315	1
Euro OTC	Tetracycline hydrochl...	120403902	1
Euro OTC	Tetracycline hydrochl...	1302012-01	10
Euro OTC	Tetracycline hydrochl...	1302012-02	17
Euro OTC	Thiamine chloride hyd...	1211008-01	1
Euro OTC	Thiamine chloride hyd...	161200401	1
Euro OTC	Thiamine chloride hyd...	1705034-01	2
Euro OTC	Thiamine chloride hyd...	1712013-01	1
Euro OTC	Thiamine chloride hyd...	1807001-01	1
Euro OTC	Thyme, rubbed	4622A130109	74
Euro OTC	Thymol	1710037-01	2
Euro OTC	Thymol	1808014-01	2
Euro OTC	Thymol	15100004-03	1
Euro OTC	Thymol	11358603	1
Euro OTC	Thymol	1510004-02	8
Euro OTC	Thymol	1410007-01	4
Euro OTC	Thymol	1510004-01	3
Euro OTC	Thymol	20151012-1	1
Euro OTC	Thymol	1202028-01	1
Euro OTC	Thymol	1510004-03	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Tobramycin sulfate	1505015-01	1
Euro OTC	Tobramycin sulfate	15I17-B04-315697	1
Euro OTC	Trace elements JK (H1...	15000526/0	1
Euro OTC	Tretinoin	1510007-01	11
Euro OTC	Tretinoin	151000	1
Euro OTC	Tretinoin	14101602	1
Euro OTC	Tretinoin	141107602	1
Euro OTC	Tretinoin	141101602	1
Euro OTC	Tretinoin	1402007-02	7
Euro OTC	Tretinoin	1411016-01	2
Euro OTC	Tretinoin	1402007-03	1
Euro OTC	Tretinoin	1402007-02	2
Euro OTC	Tretinoin	140200703	1
Euro OTC	Tretinoin	1402007	1
Euro OTC	Tretinoin	14020007-01	1
Euro OTC	Tretinoin	140007-02	1
Euro OTC	Tretinoin	140200701	1
Euro OTC	Tretinoin	130102802	1
Euro OTC	Tretinoin	30042014A	1
Euro OTC	Tretinoin	130402802	2
Euro OTC	Tretinoin	1304028-01	6
Euro OTC	Tretinoin	13041028-01	1
Euro OTC	Tretinoin	10429AX21	1
Euro OTC	Tretinoin	1008021-02	1
Euro OTC	Tretinoin	1109014-02	4
Euro OTC	Tretinoin	1402007-01	4
Euro OTC	Tretinoin	1510007-02	15
Euro OTC	Tretinoin	12059AX20	1
Euro OTC	Tretinoin	162010-02	1
Euro OTC	Tretinoin	1304028-01	1
Euro OTC	Tretinoin	1804005-01	24
Euro OTC	Tretinoin	18062082	1
Euro OTC	Tretinoin	161201001	1
Euro OTC	Tretinoin	18167302	2
Euro OTC	Tretinoin	1901001-01	6
Euro OTC	Tretinoin	190100101	2
Euro OTC	Tretinoin	161210-01	1
Euro OTC	Tretinoin	1612010-01	2
Euro OTC	Tretinoin	16004059-02	1
Euro OTC	Tretinoin	180400501	1
Euro OTC	Tretinoin	18167405	1
Euro OTC	Tretinoin	12959AX20	1
Euro OTC	Tretinoin	1604059-02	16
Euro OTC	Tretinoin	160405902	2
Euro OTC	Tretinoin	165070002	1
Euro OTC	Tretinoin	160405901	2
Euro OTC	Tretinoin	12959AX20	2
Euro OTC	Tretinoin	151000702	1
Euro OTC	Tretinoin	1604059-01	26
Euro OTC	Tretinoin	190100102	1
Euro OTC	Triamcinolone acetonide	1403009-02	13
Euro OTC	Triamcinolone acetonide	14030009-01	2
Euro OTC	Triamcinolone acetonide	141201302	1
Euro OTC	Triamcinolone acetonide	1504007-01	9
Euro OTC	Triamcinolone acetonide	15042014A	1
Euro OTC	Triamcinolone acetonide	141213-02	1
Euro OTC	Triamcinolone acetonide	1605010-01	39
Euro OTC	Triamcinolone acetonide	161006T	1
Euro OTC	Triamcinolone acetonide	1608005-01	39
Euro OTC	Triamcinolone acetonide	16D14-B11-322136	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Triamcinolone acetonide	1307043-02	2
Euro OTC	Triamcinolone acetonide	141201302	3
Euro OTC	Triamcinolone acetonide	141201301	1
Euro OTC	Triamcinolone acetonide	1407044-01	12
Euro OTC	Triamcinolone acetonide	1407044-02	13
Euro OTC	Triamcinolone acetonide	1407044-02	1
Euro OTC	Triamcinolone acetonide	3022015A	1
Euro OTC	Triamcinolone acetonide	1412013-01	12
Euro OTC	Triamcinolone acetonide	80616-5	1
Euro OTC	Triamcinolone acetonide	1412013-02	5
Euro OTC	Triamcinolone acetonide	15404007-01	1
Euro OTC	Triamcinolone acetonide	1412013-02	1
Euro OTC	Triamcinolone acetonide	NV0635M	1
Euro OTC	Triamcinolone acetonide	150800802	2
Euro OTC	Triamcinolone acetonide	1412013-01	1
Euro OTC	Triamcinolone acetonide	1601021-01	25
Euro OTC	Triamcinolone acetonide	150400701	1
Euro OTC	Triamcinolone acetonide	1504007-02	15
Euro OTC	Triamcinolone acetonide	1403009-01	4
Euro OTC	Triamcinolone acetonide	150400701	1
Euro OTC	Triamcinolone acetonide	1504007-01	1
Euro OTC	Triamcinolone acetonide	14120130	1
Euro OTC	Triamcinolone acetonide	14030009-02	1
Euro OTC	Triamcinolone acetonide	1403009-02	1
Euro OTC	Triamcinolone acetonide	1209040-01	4
Euro OTC	Triamcinolone acetonide	B0544/111404	1
Euro OTC	Triamcinolone acetonide	20130912-5	1
Euro OTC	Triamcinolone acetonide	1801020-1	1
Euro OTC	Triamcinolone acetonide	1706013	1
Euro OTC	Triamcinolone acetonide	1211013-01	8
Euro OTC	Triamcinolone acetonide	1703037-01	35
Euro OTC	Triamcinolone acetonide	1013I-02516	1
Euro OTC	Triamcinolone acetonide	1612014-01,170724T	1
Euro OTC	Triamcinolone acetonide	111005-01	1
Euro OTC	Triamcinolone acetonide	1211013-02	10
Euro OTC	Triamcinolone acetonide	12G6008-01	1
Euro OTC	Triamcinolone acetonide	1703037	3
Euro OTC	Triamcinolone acetonide	121113-02	1
Euro OTC	Triamcinolone acetonide	1703037-01	1
Euro OTC	Triamcinolone acetonide	1303001-02	8
Euro OTC	Triamcinolone acetonide	13030001-01	1
Euro OTC	Triamcinolone acetonide	170601301	2
Euro OTC	Triamcinolone acetonide	171212T/1703037-01	1
Euro OTC	Triamcinolone acetonide	1211013-02	1
Euro OTC	Triamcinolone acetonide	1307043-01	8
Euro OTC	Triamcinolone acetonide	1303001-01	2
Euro OTC	Triamcinolone acetonide	171613-01	1
Euro OTC	Triamcinolone acetonide	1612014-014	1
Euro OTC	Triamcinolone acetonide	171002-01	2
Euro OTC	Triamcinolone acetonide	1710002-1	1
Euro OTC	Triamcinolone acetonide	17100002-01	1
Euro OTC	Triamcinolone acetonide	171302-01	1
Euro OTC	Triamcinolone acetonide	180102001	1
Euro OTC	Triamcinolone acetonide	1710002-01,180626T	1
Euro OTC	Triamcinolone acetonide	1706013-01	45
Euro OTC	Triamcinolone acetonide	181016181	1
Euro OTC	Triamcinolone acetonide	PT49	1
Euro OTC	Triamcinolone acetonide	1209040-02	3
Euro OTC	Triamcinolone acetonide	1307043-01	6
Euro OTC	Triamcinolone acetonide	180500101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Triamcinolone acetonide	1602005-01	1
Euro OTC	Triamcinolone acetonide	1608005	1
Euro OTC	Triamcinolone acetonide	1612014-01	56
Euro OTC	Triamcinolone acetonide	1307043-02	13
Euro OTC	Triamcinolone acetonide	1206008-01	5
Euro OTC	Triamcinolone acetonide	180500102	1
Euro OTC	Triamcinolone acetonide	1903015-01	2
Euro OTC	Triamcinolone acetonide	181001501	1
Euro OTC	Triamcinolone acetonide	1805001	1
Euro OTC	Triamcinolone acetonide	161214-01	1
Euro OTC	Triamcinolone acetonide	1805001-01	12
Euro OTC	Triamcinolone acetonide	1110005-01	2
Euro OTC	Triamcinolone acetonide	1112023-01	1
Euro OTC	Triamcinolone acetonide	1112023-02	4
Euro OTC	Tricalcium phosphate	1206002-01	3
Euro OTC	Tricalcium phosphate	1304025-01	1
Euro OTC	Triclosan	169005-02	1
Euro OTC	Triclosan	140603201	1
Euro OTC	Triclosan	170113T1609005-01	1
Euro OTC	Triclosan	1401038-01	12
Euro OTC	Triclosan	1701032-02	16
Euro OTC	Triclosan	190100601	1
Euro OTC	Triclosan	170511-01	1
Euro OTC	Triclosan	1505009-01	14
Euro OTC	Triclosan	CHB1303011-01	1
Euro OTC	Triclosan	15058946	2
Euro OTC	Triclosan	1307050-01	8
Euro OTC	Triclosan	1550229	1
Euro OTC	Triclosan	16047002-01	1
Euro OTC	Triclosan	5860088	1
Euro OTC	Triclosan	180149-01	1
Euro OTC	Triclosan	1245	1
Euro OTC	Triclosan	1501032-01	10
Euro OTC	Triclosan	130301101	1
Euro OTC	Triclosan	190100602	1
Euro OTC	Triclosan	170825-01	1
Euro OTC	Triclosan	1406032-01	12
Euro OTC	Triclosan	171121T/1705011-01	1
Euro OTC	Triclosan	1701032_01	1
Euro OTC	Triclosan	1708025_01	1
Euro OTC	Triclosan	1609005-02	22
Euro OTC	Triclosan	1408024-01	12
Euro OTC	Triclosan	15150229	2
Euro OTC	Triclosan	1801049	1
Euro OTC	Triclosan	120603402	1
Euro OTC	Triclosan	12276024	1
Euro OTC	Triclosan	170132-01	1
Euro OTC	Triclosan	1901006-01	6
Euro OTC	Triclosan	1708025-01	25
Euro OTC	Triclosan	140103801	1
Euro OTC	Triclosan	1609005-01	25
Euro OTC	Triclosan	1507041-01	8
Euro OTC	Triclosan	170103202	1
Euro OTC	Triclosan	1202018-02	2
Euro OTC	Triclosan	170526T/1701032-02	1
Euro OTC	Triclosan	1211025-01	11
Euro OTC	Triclosan	120203-02	1
Euro OTC	Triclosan	1510045-01	15
Euro OTC	Triclosan	150500901	1
Euro OTC	Triclosan	1801049-21	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Triclosan	17080-01	1
Euro OTC	Triclosan	1701032	1
Euro OTC	Triclosan	160700201	1
Euro OTC	Triclosan	700228	1
Euro OTC	Triclosan	1607002-01	36
Euro OTC	Triclosan	160702-01	1
Euro OTC	Triclosan	160531T	1
Euro OTC	Triclosan	1501032-02	6
Euro OTC	Triclosan	4486611	2
Euro OTC	Triclosan	160200301	1
Euro OTC	Triclosan	1206034-01	1
Euro OTC	Triclosan	160900502	2
Euro OTC	Triclosan	1602003	1
Euro OTC	Triclosan	1206034-02	5
Euro OTC	Triclosan	1602003-01	29
Euro OTC	Triclosan	1303011-01	5
Euro OTC	Triclosan	12276023	1
Euro OTC	Trometamol	1203005-02	3
Euro OTC	Tyrosine	L1210050-01	1
Euro OTC	Tyrosine	L1401037	1
Euro OTC	Tyrosine	L1506022-02	1
Euro OTC	Ubidecarenone	L1403020-01	1
Euro OTC	Ubidecarenone	1403020-02	1
Euro OTC	Ubidecarenone	1902030-01	1
Euro OTC	Ubidecarenone	1515M-07029	1
Euro OTC	Ubidecarenone	1014E-07029	1
Euro OTC	Ubidecarenone	1605030-01	3
Euro OTC	Ubidecarenone	1802025-01	2
Euro OTC	Urea	1212007-01Mg:300g	1
Euro OTC	Urea	16I13-B02-332953	2
Euro OTC	Urea	UP250171017	1
Euro OTC	Urea	180400701	1
Euro OTC	Urea	1212007-01	8
Euro OTC	Urea	1804007-01	19
Euro OTC	Urea	1404024-01	14
Euro OTC	Urea	1212007-02	3
Euro OTC	Urea	170701501	1
Euro OTC	Urea	1707015-01	25
Euro OTC	Urea	72/14	1
Euro OTC	Urea	V076-2016	1
Euro OTC	Urea	1110023-01	2
Euro OTC	Urea	13064239	3
Euro OTC	Urea	1503042-01	25
Euro OTC	Urea	1811012-01	6
Euro OTC	Urea	150304201	2
Euro OTC	Urea	1515404	2
Euro OTC	Urea	18040701	1
Euro OTC	Urea	1601044-01	27
Euro OTC	Urea	1604044-01	26
Euro OTC	Urea	160404401	2
Euro OTC	Urea	1702023-01	24
Euro OTC	Urea	1804007-02	2
Euro OTC	Urea	170202301	3
Euro OTC	Urea	14107101	12
Euro OTC	Urea	160104401	1
Euro OTC	Urea	1804007-1	1
Euro OTC	Valine	1710003-01	1
Euro OTC	Valine	1604055-01	2
Euro OTC	Valine	L1309029-02	1
Euro OTC	Vitamin B6 (H1103)	15000279/1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Euro OTC	Wax, bleached	16163104	1
Euro OTC	Wax, yellow	1503013-01	1
Euro OTC	Wax, yellow	1407013-01	2
Euro OTC	Wax, yellow	1201012-01	1
Euro OTC	Wax, yellow	18113505	13
Euro OTC	Xuan Shen	H0537021	1
Euro OTC	Xylose	16206405-	1
Euro OTC	Zinc gluconate	27061711	1
Euro OTC	Zinc gluconate	1606018-02	36
Euro OTC	Zinc gluconate	1807019-01	14
Euro OTC	Zinc gluconate	1807019	1
Euro OTC	Zinc gluconate	1709012-0103707875	1
Euro OTC	Zinc gluconate	1701017-01	31
Euro OTC	Zinc gluconate	1312010-01	1
Euro OTC	Zinc gluconate	1706009-01	66
Euro OTC	Zinc gluconate	170901201	9
Euro OTC	Zinc gluconate	1711004-01	27
Euro OTC	Zinc gluconate	1405011-01	1
Euro OTC	Zinc gluconate	1606018-01	4
Euro OTC	Zinc sulfate heptahyd...	1201031-02	1
Euro OTC	Zinc sulfate heptahyd...	1201031-2	1
Euro OTC	Zinc sulfate heptahyd...	1309011_0	1
Euro OTC	Zinc sulfate heptahyd...	1105002-01	1
Euro OTC	Zinc sulfate heptahyd...	1201031-01	1
Euro OTG	Pilocarpine hydrochlo...	1305044-01	2
Euro OTG/ AHD	Tetracaine hydrochlor...	1706018-01	10
euro/ Anzag	Triamcinolone acetone...	3315I02516	1
Euro/ Phönix	Estradiol hemihydrate	3417M-02845	1
euro/ Sanacorp	Camphor	17050602	2
Euro/ Sanacorp	Metronidazole	150800402	1
Euro/ Sanacorp	Progesterone, microni...	C01-131116(M)	1
euro/ Noweda	Betamethasone dipropi...	1705030-01	1
EUro/AHD	Aluminium chloride he...	121001401	3
EUro/AHD	Caffeine	121012-01	1
EUro/AHD	Caffeine	1307049-01	3
EUro/AHD	Calcium chloride dihy...	1205014-02	2
EUro/AHD	Chlortetracycline hyd...	1108051-02	11
EUro/AHD	Prednisolone acetate	1306040-02	16
Euro/Alliance Healthc...	Lidocaine hydrochlori...	150523-01	1
Euro/Alliance Healthc...	Sodium tetraborate	171100201	2
Euro/Anzag	Minoxidil	150903501	2
Euro/Gehe	Emulsifying cetostear...	150401-01	1
Euro/Gehe	Fructose	1604024-01	6
Euro/Gehe	Fructose	1607030-01	17
Euro/Gehe	Minoxidil	140528-01	2
Euro/Noweda	Ambroxol hydrochloride	1509032-01	2
Euro/Noweda	Cellulose, microcryst...	1410002-02	1
Euro/Noweda	Triamcinolone acetone...	161201401	1
Euro/Phönix	Carbomers 35000 / 50000	161002401	1
Euro/Phönix	Estriol	1509042-01	1
Euro/Phönix	Lidocaine hydrochlori...	141101010-01	1
Euro/Phönix	Minoxidil	161201101	3
Euro/Phönix	Norethisterone acetate	1306033-01	1
Euro/Phönix	Triclosan	170103201	1
EuRo/Phönix	Chlortetracycline hyd...	160401501	1
Euro/Sanacorp	alpha-tocopheryl acet...	150700101	1
Euro/Sanacorp	Carbomers 35000 / 50000	1305007-01	2
Euro/Sanacorp	Metronidazole	1609028	3
Euro/Sanacorp	Wax, yellow	1708506	1
euroct/gehe	L-tryptophan	1212010-01	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
euroct/gehe	Triamcinolone acetonide	1403009	1
euroct/Gehe	Prednisolone, microni...	1212021-01	19
euroho/gehe	Sodium hydrogen carbo...	123030-02	1
Euroo	Glucose monohydrate	1604006-01	25
Euroo/Anzag	Calcium pantothenate	1509024-01	6
Euroo/direkt	Triclosan	1705011-01	38
Euroo/Phoenix	Aluminium chloride he...	120014-02	1
Euroo/Phoenix	Povidone iodine	G3115BPZ0	1
Euroo/Sana	Povidone iodine	1501029-01	3
Europ Ots Pharma	Betamethasone valerate	12111-N01	1
Europ Ots Pharma	Triclosan	13D02-N02	2
Europharm	Menthol	130501-01	1
Europharm/Otto	Diclofenac sodium	135796	1
Europharm/Otto	Lactose monohydrate	1106022-01	1
Europharm/Otto	Prednisolone, microni...	1108055-04	21
Europharm/Otto	Tetracycline hydrochl...	1204039-02	16
Europharm/Sanacorp	Triamcinolone acetonide	B0544/111040	2
europharma	Prednisolone, microni...	1002001240	1
EuroPharma	Metronidazole	141201601	1
EuroPTC/Phönix	Miconazole nitrate	14149316	1
EuroTOC	Ambroxol hydrochloride	1602015-01	2
EuroTOC / Noweda	Sodium benzoate	130902701	1
EuroTOC / Noweda	Tetracaine hydrochlor...	131000601	2
EurRho/ Anzag	Diclofenac sodium	1403013-03	5
f	Hydrocortisone acetate	333961	1
F	5-aminolevulinic acid...	13C-19-N07	1
fa	Citric acid	13K18-B34-291523	5
Fa Roth	Histamine dihydrochlo...	47252901	1
Facron	Macrogols 4000	7022014C	1
Fag / Noweda	Dequalinium chloride	15B16-B04-321326	2
Fagaron	Salicylic acid	14Lo2-B02-303558	1
Fagon	Mannitol	13K20-B01-292064	1
Fagraon	Bacitracin	14E23-B03-303424	3
Fagrin/Noweda	Ascorbyl palmitate	18B12-B09525720	1
Fagro	Carmellose sodium (40...	17C20-B07-337104	1
Fagro	Carmellose sodium (40...	17K06-B03-345975	1
Fagro	Carmellose sodium (40...	16G05-B05-329314	1
Fagro	Carmellose sodium (40...	15G03-B02-321790	1
Fagro	Emulsifying cetostear...	5211Q-01078	1
Fagro	Metoprolol tartrate	16MD000095	1
Fagro	Mometasone furoate, m...	1511B02	1
Fagro WE:19.05.16 EK:...	L-isoleucine	15A14-B09-314881	1
Fagro/Anzag	Oxytetracycline hydro...	14122-B01-301464	1
Fagro/Gehe	Hydrocortisone acetate	15A07-B03_306545	1
Fagro/Gehe	Lidocaine	15a15b01304284	1
Fagro/Sanacorp	Oxybutynin hydrochlor...	16/23-B03-32891613	1
Fagron	(R,S)-epinephrine hyd...	16H23-B98-325052	1
Fagron	(R,S)-epinephrine hyd...	10022015E	1
Fagron	(R,S)-epinephrine hyd...	2021503	1
Fagron	(R,S)-epinephrine hyd...	14A22-B05	1
Fagron	5-aminolevulinic acid...	14K13-B01	1
Fagron	5-aminolevulinic acid...	14K04-B01-297971	2
Fagron	5-aminolevulinic acid...	13H13-N07	2
Fagron	5-aminolevulinic acid...	11I30-N01	1
Fagron	5-aminolevulinic acid...	7447	1
Fagron	5-aminolevulinic acid...	15C06-B02-305245	1
Fagron	5-aminolevulinic acid...	13C19-N05	3
Fagron	5-aminolevulinic acid...	12K26-N05	1
Fagron	5-aminolevulinic acid...	12H20-N01	1
Fagron	5-aminolevulinic acid...	12c23-n05	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	5-aminolevulinic acid...	13J29-N03	2
Fagron	5-aminolevulinic acid...	17F01-B06-337164	2
Fagron	5-aminolevulinic acid...	16A11-B02-318406	1
Fagron	5-aminolevulinic acid...	17F01-B05-338204	1
Fagron	5-aminolevulinic acid...	16K17-B09-328197	1
Fagron	5-aminolevulinic acid...	1701-b05-342090	1
Fagron	5-aminolevulinic acid...	17A12-B07-331744	1
Fagron	5-aminolevulinic acid...	17122-B02-343410	1
Fagron	5-aminolevulinic acid...	17A12-B07-331745	1
Fagron	5-aminolevulinic acid...	17/22-B02-342555	1
Fagron	5-aminolevulinic acid...	16G27-B10	1
Fagron	5-aminolevulinic acid...	15J23-B06	1
Fagron	5-aminolevulinic acid...	17L18-B05-344673	1
Fagron	5-aminolevulinic acid...	18C16-B05-347714	1
Fagron	5-aminolevulinic acid...	18G13-B06	1
Fagron	5-aminolevulinic acid...	16A11-B02	1
Fagron	5-aminolevulinic acid...	15J23-B06-314284	1
Fagron	5-aminolevulinic acid...	16I12-B06	1
Fagron	5-aminolevulinic acid...	16A11-B02-316734	1
Fagron	5-aminolevulinic acid...	16A11-B02-317459	3
Fagron	5-aminolevulinic acid...	18C16-B05	1
Fagron	5-aminolevulinic acid...	17F01-B06_337165	1
Fagron	5-aminolevulinic acid...	13K27-N02	3
Fagron	8-quinolinol sulfate ...	13G10-N06	1
Fagron	8-quinolinol sulfate ...	11L08-N07	2
Fagron	8-quinolinol sulfate ...	10F01-N15	1
Fagron	8-quinolinol sulfate ...	13/21-N23	1
Fagron	8-quinolinol sulfate ...	12A04-N05	1
Fagron	8-quinolinol sulfate ...	12D26-N04	1
Fagron	8-quinolinol sulfate ...	12H21-N08	1
Fagron	8-quinolinol sulfate ...	11I09-N09	1
Fagron	8-quinolinol sulfate ...	11/09-N09	1
Fagron	8-quinolinol sulfate ...	13032014C	1
Fagron	8-quinolinol sulfate ...	13D24-N11	1
Fagron	8-quinolinol sulfate ...	13A04-N05	1
Fagron	8-quinolinol sulfate ...	11L19-N14	1
Fagron	Acetylcysteine	13A22N07	1
Fagron	Acetylsalicylic acid	14B10-B02-296841	1
Fagron	Aciclovir	17G11-B01-342081	1
Fagron	Aciclovir	15AS15-B10-307743	1
Fagron	Aciclovir	15A15-B10-307743	1
Fagron	Aciclovir	146501524	1
Fagron	Aciclovir	16K09-B05-328196	4
Fagron	Aciclovir	18I18-B04-355420	1
Fagron	Aciclovir	15A15-B10-315379	13
Fagron	Aciclovir	17G11-B01-339460	1
Fagron	Aciclovir	15A15-B10-323362	1
Fagron	Aciclovir	14g07-b05-296462	1
Fagron	Albendazole	14B05-B41-305249	1
Fagron	Albendazole	14b05-b41-314321	1
Fagron	Albendazole	14B05-B41-314321	1
Fagron	Albendazole	14B05-B41	1
Fagron	Albendazole	14B05B41-312959	1
Fagron	Alchemilla, cut	11182802	1
Fagron	Alfatradiol	13L30-B32	1
Fagron	Alfatradiol	14/22-B30-302017	1
Fagron	Alfatradiol	13G19-N06	2
Fagron	Alfatradiol	15G07-B01	1
Fagron	Alfatradiol	15J27-B08-321786	1
Fagron	Alfatradiol	15G07-B01-310228	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Alfatradiol	15J27-B08	1
Fagron	Alfatradiol	14G07-B99-294945	3
Fagron	Alfatradiol	14G07-B99	1
Fagron	Alfatradiol	13L30B32	1
Fagron	Alfatradiol	13L30-B32-290739	1
Fagron	Alfatradiol	13L30-B32	1
Fagron	Alfatradiol	13L30-B32- . . .	1
Fagron	Alfatradiol	13130b32290742	1
Fagron	Alfatradiol	15J27-B08-318783	1
Fagron	Alfatradiol	13L30-B32-290742	1
Fagron	Alfatradiol	14722-B30-302018	1
Fagron	Alfatradiol	5102432	1
Fagron	Alfatradiol	13G19-N09	1
Fagron	Alfatradiol	13G19-No6	1
Fagron	Alfatradiol	14I22-B30	1
Fagron	Alfatradiol	15G07-B01-310227	5
Fagron	Alfatradiol	15G07-B01-311881	2
Fagron	Alfatradiol	16g05-b01-324614	1
Fagron	Alfatradiol	13G19-N06	1
Fagron	Alfatradiol	324614	1
Fagron	Alfatradiol	16G05-Bo1-332379	1
Fagron	Alfatradiol	16G05-B01-334055	1
Fagron	Alfatradiol	17G03-B11	2
Fagron	Alfatradiol	16G05-B01-337352	1
Fagron	Alfatradiol	17G03-B11-345093	2
Fagron	Alfatradiol	11B07-N05	1
Fagron	Alfatradiol	17G03-B11-346962	1
Fagron	Alfatradiol	18C23-B05	2
Fagron	Alfatradiol	9011904	1
Fagron	Alfatradiol	18C23-B05-352724	1
Fagron	Alfatradiol	18C23-B05-354802	1
Fagron	Alfatradiol	12D19-N03	6
Fagron	Alfatradiol	14G07-B99-294945	1
Fagron	Alfatradiol	17G03-B11-348863	1
Fagron	Alfatradiol	11K25-N01	3
Fagron	Alfatradiol	11D07-N19	1
Fagron	Alfatradiol	11G22-N11	2
Fagron	Alfatradiol	12J27-B08-314546	1
Fagron	Alfatradiol	13L30B32-290741	1
Fagron	Alfatradiol	16G05-801	1
Fagron	Alfatradiol	16G05-B01-325108	1
Fagron	Alfatradiol	14G07-B99-294942	1
Fagron	Alfatradiol	12L19-No8	1
Fagron	Alfatradiol	11k25-N01	1
Fagron	Alfatradiol	16G05-B01-328254	2
Fagron	Alfatradiol	13G19-No6	1
Fagron	Alfatradiol	14C20-B02-294943	4
Fagron	Alfatradiol	13G19-N06	18
Fagron	Alfatradiol	16G05-B01-326866	1
Fagron	Alfatradiol	16G05-B01-335789	5
Fagron	Alfatradiol	XP101UT	1
Fagron	Alfatradiol	16G05-B01-324614	2
Fagron	Alfatradiol	15J27B08318164	1
Fagron	Alfatradiol	15J27-B08-314546	1
Fagron	Alfatradiol	15G07-B01-310229	3
Fagron	Alfatradiol	14C20-B02-294943	1
Fagron	Alfatradiol	15J27-B08-318164	3
Fagron	Allantoin	13A21-N07	3
Fagron	Allantoin	13I03-B05-289430	1
Fagron	Allantoin	11J11-N11	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Allantoin	12B09-N09	1
Fagron	Allantoin	13G16-N07	7
Fagron	Allantoin	5011E-01087	1
Fagron	Allantoin	12C30-N04	1
Fagron	Allopurinol	14E15-B08	1
Fagron	Allopurinol	14E15-B05-315557	13
Fagron	Allopurinol	14E15-B05-320787	1
Fagron	Almonds, bitter, whole	40A161005	4
Fagron	alpha-tocopheryl acet...	13H30-B02-287169	1
Fagron	alpha-tocopheryl acet...	16H30-B05-326792	1
Fagron	Aluminium chloride he...	11G05-No2	1
Fagron	Aluminium chloride he...	11I23-N06	3
Fagron	Aluminium chloride he...	12I27-N01	4
Fagron	Aluminium chloride he...	2913U1095	1
Fagron	Aluminium chloride he...	13J31-B06-287615	1
Fagron	Aluminium chloride he...	12082013D	1
Fagron	Aluminium chloride he...	20122013E	1
Fagron	Aluminium chloride he...	13I26-N13	2
Fagron	Aluminium chloride he...	13I26-N13	1
Fagron	Aluminium chloride he...	12047512	3
Fagron	Aluminium chloride he...	5213A1095	1
Fagron	Aluminium chloride he...	105282-0001	1
Fagron	Aluminium chloride he...	13d22no9	1
Fagron	Aluminium chloride he...	11G05-N02	1
Fagron	Aluminium chloride he...	11/23-N06	2
Fagron	Aluminium chloride he...	11 23-N06	1
Fagron	Aluminium chloride he...	C012A	1
Fagron	Aluminium chloride he...	12/27-N01	7
Fagron	Aluminium chloride he...	C181A	2
Fagron	Aluminium chloride he...	12I27-N01	11
Fagron	Aluminium chloride he...	1910211095	1
Fagron	Aluminium chloride he...	12G05-N10	23
Fagron	Aluminium chloride he...	12/11-No3	1
Fagron	Aluminium chloride he...	11/23-N06B206A	1
Fagron	Aluminium chloride he...	12L11-N03	3
Fagron	Aluminium chloride he...	12C07-N01	7
Fagron	Aluminium chloride he...	8547287	1
Fagron	Aluminium chloride he...	2013U-01095	1
Fagron	Aluminium chloride he...	12i27-N01	2
Fagron	Aluminium chloride he...	12/27-n01	1
Fagron	Aluminium chloride he...	13D22N09	2
Fagron	Aluminium chloride he...	13D22-N09	26
Fagron	Aluminium chloride he...	C145A	1
Fagron	Aluminium chloride he...	12/27/N01	1
Fagron	Aluminium chloride he...	1207211095	1
Fagron	Aluminium potassium s...	103638	1
Fagron	Aluminium potassium s...	12/24-M02	1
Fagron	Aluminium potassium s...	A1-063-AL	1
Fagron	Aluminium potassium s...	13D29-N04	1
Fagron	Aluminium potassium s...	12I24-M03	1
Fagron	Aluminium potassium s...	A1-063-A1	1
Fagron	Aluminium potassium s...	13D08-N01	1
Fagron	Ambroxol hydrochloride	12K12-N02	1
Fagron	Ambroxol hydrochloride	16G12-B06-326184	1
Fagron	Amifampridine	13L12-M03	1
Fagron	Amifampridine	16D05-B01-318216	1
Fagron	Amifampridine	11/14-N03	1
Fagron	Amifampridine	14H11-B06	3
Fagron	Amifampridine	16D05-B01	1
Fagron	Amifampridine	11K07-N02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Amifampridine	15B20-B02	2
Fagron	Amifampridine	13K14-M01	1
Fagron	Amifampridine	16K15-B03-330546	1
Fagron	Amifampridine	13F07-N07	1
Fagron	Amifampridine	13L12-M06	2
Fagron	Amiloride hydrochloride	17J02-B04	1
Fagron	Amiloride hydrochloride	17J02-B04-341175	1
Fagron	Amitriptyline hydroch...	19C08-B02-362281	1
Fagron	Amitriptyline hydroch...	17E11-B02-340187	1
Fagron	Amitriptyline hydroch...	19C08-B02-362280	1
Fagron	Amitriptyline hydroch...	18A16-B06-348643	1
Fagron	Amitriptyline hydroch...	18A16-B06-344853	2
Fagron	Amitriptyline hydroch...	18A16-B06-353065	2
Fagron	Amphotericin B	5173	1
Fagron	Amphotericin B	16D05-B04	1
Fagron	Aniseed	11370304	6
Fagron	Arginine	16H24-B04-335673	2
Fagron	Arginine	17G10-B07	1
Fagron	Arginine	13F14-N02	2
Fagron	Arginine	12C21-N04	1
Fagron	Arginine	16032361	1
Fagron	Arginine	14/24-B08-303720	4
Fagron	Arginine	13031801	1
Fagron	Arginine	16H24-B04-334063	1
Fagron	Arginine	16D11-B01-323566	1
Fagron	Arginine	16H24-B04-335674	7
Fagron	Arginine	14062019	1
Fagron	Arginine	14/24-B08-303901	1
Fagron	Arginine	13J28-B02-290463	7
Fagron	Arginine	14092035	2
Fagron	Arginine	1424-B08-303720	1
Fagron	Arginine	14I24-B08-303901	1
Fagron	Arginine	15C10-B03-318784	1
Fagron	Arginine	16H24-B04-328914	1
Fagron	Arginine	16D11-B01-322156	1
Fagron	Arginine	15102258	1
Fagron	Arginine	14I24-B-303720	1
Fagron	Arginine	14B25-B04	2
Fagron	Arginine (L-) (H1036)	1324044301	2
Fagron	Arginine hydrochloride	13J28-B03	2
Fagron	Arginine hydrochloride	13j28-B03-287768	1
Fagron	Arginine hydrochloride	13j28-b03-287768	1
Fagron	Arginine hydrochloride	13J28-B032x	1
Fagron	Arginine hydrochloride	13J28-B03-287766	1
Fagron	Arginine hydrochloride	12C20-N08	1
Fagron	Ascorbyl palmitate	17H04-B01-340474	2
Fagron	Ascorbyl palmitate	K50642890	1
Fagron	Ascorbyl palmitate	16 30-B03-525358	1
Fagron	Ascorbyl palmitate	16/30-B03-330625	2
Fagron	Ascorbyl palmitate	16/30-B03	1
Fagron	Ascorbyl palmitate	16I30-B03-329375	1
Fagron	Ascorbyl palmitate	16I30-B03-330625	1
Fagron	Ascorbyl palmitate	16I30-B03-326543	1
Fagron	Ascorbyl palmitate	16B10-B04-525228	1
Fagron	Ascorbyl palmitate	16I30-B03-332623	3
Fagron	Ascorbyl palmitate	15F01-B99-525121	2
Fagron	Ascorbyl palmitate	16I30-B03-525358	2
Fagron	Ascorbyl palmitate	15F01-B99-525063	1
Fagron	Ascorbyl palmitate	16B10-B04-525199	3
Fagron	Ascorbyl palmitate	15F01-B99-525122	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Ascorbyl palmitate	15f01-b99-309774	1
Fagron	Ascorbyl palmitate	16I30-B03	2
Fagron	Ascorbyl palmitate	15F01-B99-52063	1
Fagron	Ascorbyl palmitate	19D25-B04	2
Fagron	Ascorbyl palmitate	16/30-B03-329375	1
Fagron	Ascorbyl palmitate	14I03-N01	1
Fagron	Ascorbyl palmitate	17E18-B06-525524	3
Fagron	Ascorbyl palmitate	17e18-b06	2
Fagron	Ascorbyl palmitate	17H04-B01-338484	7
Fagron	Ascorbyl palmitate	525426	1
Fagron	Ascorbyl palmitate	525489	1
Fagron	Ascorbyl palmitate	525525	1
Fagron	Ascorbyl palmitate	525524	1
Fagron	Ascorbyl palmitate	13C06-B40-291206	2
Fagron	Ascorbyl palmitate	17H04-B01-525632	1
Fagron	Ascorbyl palmitate	16,30-B03	1
Fagron	Ascorbyl palmitate	K48951290	1
Fagron	Ascorbyl palmitate	17E18-B06-525514	1
Fagron	Ascorbyl palmitate	17E18-B06-525525	1
Fagron	Ascorbyl palmitate	15F01-B99	3
Fagron	Ascorbyl palmitate	18B12-B09-525719	2
Fagron	Ascorbyl palmitate	K48951390	1
Fagron	Ascorbyl palmitate	K45532790	2
Fagron	Ascorbyl palmitate	23061810	1
Fagron	Ascorbyl palmitate	17e18b06334914	1
Fagron	Ascorbyl palmitate	18A10B03344364	1
Fagron	Ascorbyl palmitate	16I30-B03-525380	2
Fagron	Ascorbyl palmitate	18B12B09525719	1
Fagron	Ascorbyl palmitate	18J22-B03-357638	1
Fagron	Ascorbyl palmitate	18J22-B03-355582	2
Fagron	Ascorbyl palmitate	16L30-B03-332623	1
Fagron	Ascorbyl palmitate	18B12-B09-525836	1
Fagron	Ascorbyl palmitate	16I30-B03-525426	1
Fagron	Ascorbyl palmitate	10061708	1
Fagron	Ascorbyl palmitate	18A10-B03-344364	10
Fagron	Ascorbyl palmitate	925749	1
Fagron	Ascorbyl palmitate	K47747490	1
Fagron	Ascorbyl palmitate	17H04-B01	1
Fagron	Ascorbyl palmitate	1630-B03-525489	1
Fagron	Ascorbyl palmitate	17E18-B06	2
Fagron	Ascorbyl palmitate	18B12-B09-348867	5
Fagron	Bacitracin	13C13-N10	1
Fagron	Bacitracin	211531	1
Fagron	Bacitracin	13B25-B01	1
Fagron	Bacitracin	13A03-N11	1
Fagron	Bacitracin	14A21-B09	1
Fagron	Bacitracin	14A21-B09-290645	2
Fagron	Bacitracin	12G12-N17	1
Fagron	Bacitracin	13D26-N08	1
Fagron	Bacitracin	14A21-B09-290643	2
Fagron	Bacitracin	17L11-B02-343686	1
Fagron	Bacitracin	11J17-N01	2
Fagron	Bacitracin	17L11-B02	3
Fagron	Bacitracin	16B17-B03	2
Fagron	Bacitracin	17L11-B02-343685	4
Fagron	Bacitracin	1612-B033-334052	1
Fagron	Bacitracin	14A21-B09-292908	1
Fagron	Bacitracin	16B17-B03-322939	1
Fagron	Bacitracin	16B17-B03-322891	5
Fagron	Bacitracin	16/12-B03-329401	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Bacitracin	16/12B03335331	1
Fagron	Bacitracin	16b17-b03-322939	1
Fagron	Bacitracin	16I12-B03-334052	2
Fagron	Bacitracin	16/12-B03-331152	1
Fagron	Bacitracin	2917M-03314	1
Fagron	Bacitracin	16I12-B03-335331	1
Fagron	Bacitracin	16I12-B03-329401	1
Fagron	Bacitracin	17L11-802	1
Fagron	Bacitracin	15C12-B06	1
Fagron	Bacitracin	17L11-B02-346531	1
Fagron	Bacitracin	C2280017	3
Fagron	Bacitracin	14E23-B03-296243	1
Fagron	Bacitracin	17L11-B02-35400	1
Fagron	Bacitracin	14E23-B03-296243	7
Fagron	Bacitracin	14A21-B09-290643	1
Fagron	Bacitracin	14E23-B03-296245	2
Fagron	Bacitracin	14A21-B09	1
Fagron	Bacitracin	16I12-B03,170913B	1
Fagron	Bacitracin	14E23-B03-333424	1
Fagron	Bacitracin	14E23-B03	3
Fagron	Bacitracin	17L11B02357635	1
Fagron	Bacitracin	14e23-b03-296243	1
Fagron	Bacitracin	14E23B03-303424	1
Fagron	Bacitracin	17L11B02348461	1
Fagron	Bacitracin	16I12-B03	2
Fagron	Beclometasone dipropi...	13E02-B40-293328	2
Fagron	Beclometasone dipropi...	13I05-N03	3
Fagron	Beclometasone dipropi...	18C05-B01-347367	1
Fagron	Beclometasone dipropi...	12D12-N03	4
Fagron	Beclometasone dipropi...	15F29-B09-314376	1
Fagron	Beclometasone dipropi...	16B16-B03-322517	1
Fagron	Beclometasone dipropi...	16B16-B03-318877	1
Fagron	Beclometasone dipropi...	13/05-N03	1
Fagron	Beclometasone dipropi...	14D04-B04-302335	2
Fagron	Beclometasone dipropi...	14D04-B04-302336	2
Fagron	Beclometasone dipropi...	14D04-B04-302336	1
Fagron	Beclometasone dipropi...	13/05-N03	1
Fagron	Beclometasone dipropi...	11L01-N03	1
Fagron	Beclometasone dipropi...	13E02-B40-293327	3
Fagron	Beclometasone dipropi...	13I05N03	1
Fagron	Beclometasone dipropi...	13E02-B40-296354	2
Fagron	Beclometasone dipropi...	11F28-N02	2
Fagron	Beclometasone dipropi...	12G04-N11	3
Fagron	Benzocaine	11/22-N3	1
Fagron	Benzocaine	13I06-B02-300494	1
Fagron	Benzocaine	14G02-B01-310324	7
Fagron	Benzocaine	15C19-B04-316239	2
Fagron	Benzocaine	18F29B07190519	1
Fagron	Benzocaine	18F18-B50	2
Fagron	Benzocaine	581015	1
Fagron	Benzocaine	18F18-B50-352167	2
Fagron	Benzocaine	15C19-B04	2
Fagron	Benzocaine	12E014-N06	1
Fagron	Benzocaine	14G02-B01-318318	2
Fagron	Benzocaine	15C19-B04-322019	6
Fagron	Benzocaine	15C-19-B04-325359	1
Fagron	Benzocaine	13F06N10	1
Fagron	Benzocaine	14G02B01318318	1
Fagron	Benzocaine	16D13	1
Fagron	Benzocaine	16D13-B01-337264	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Benzocaine	16D13-B01-342411	3
Fagron	Benzocaine	16D13-B01	1
Fagron	Benzocaine	16D13-B01-345969	2
Fagron	Benzocaine	16D13B01341861	2
Fagron	Benzocaine	16D113-B01-345969	1
Fagron	Benzocaine	16D13-B01-348062	3
Fagron	Benzocaine	16D13-B01-334195	1
Fagron	Benzocaine	13L06Bo2	1
Fagron	Benzocaine	13L06-B02-300494	4
Fagron	Benzocaine	15C19-B04-308254	1
Fagron	Benzocaine	4413E-01118	1
Fagron	Benzocaine	10K02-N01	1
Fagron	Benzocaine	12E04-N06	9
Fagron	Benzocaine	11/22-N03	1
Fagron	Benzocaine	BZ-27/1112	1
Fagron	Benzocaine	11K08-B40	4
Fagron	Benzocaine	104210-0001	1
Fagron	Benzocaine	11k08-B40-291841	1
Fagron	Benzocaine	11K08-B40-291841	5
Fagron	Benzocaine	BZ-08/11-12	1
Fagron	Benzocaine	11K08-B40-291840	1
Fagron	Benzocaine	11k08-b40-291841	1
Fagron	Benzocaine	13L06-B02	5
Fagron	Benzocaine	18F18-B50-352171	2
Fagron	Benzocaine	13L06-B02-294100	2
Fagron	Benzocaine	13L06-B02_297505	1
Fagron	Benzocaine	400313	1
Fagron	Benzydamine hydrochlo...	15A14-B05	1
Fagron	Benzydamine hydrochlo...	15I17-B02	1
Fagron	Benzydamine hydrochlo...	15A14-B05-303068	1
Fagron	Benzydamine hydrochlo...	13A15-N03	1
Fagron	Benzydamine hydrochlo...	12C20-N02	1
Fagron	Benzydamine hydrochlo...	17H01-B03-348469	1
Fagron	Benzydamine hydrochlo...	12K07-N16	2
Fagron	Benzydamine hydrochlo...	1503309P	1
Fagron	Benzydamine hydrochlo...	17H01-B03-345051	1
Fagron	Benzydamine hydrochlo...	17H01-B03_34641	1
Fagron	Benzydamine hydrochlo...	17H01-B03-345050	3
Fagron	Benzydamine hydrochlo...	1703305P	1
Fagron	Benzydamine hydrochlo...	16D05-B50-332334	1
Fagron	Benzydamine hydrochlo...	16d05-b50-332333	1
Fagron	Benzydamine hydrochlo...	16D05-b50-332334	1
Fagron	Benzydamine hydrochlo...	23.03.2017	1
Fagron	Benzydamine hydrochlo...	3280417	1
Fagron	Benzydamine hydrochlo...	16B09-B05-318272	2
Fagron	Benzydamine hydrochlo...	16B09-B05-320007	1
Fagron	Benzydamine hydrochlo...	12J10-N03	1
Fagron	Benzydamine hydrochlo...	1503313P	1
Fagron	Benzydamine hydrochlo...	14B25-B01-294054	1
Fagron	Benzydamine hydrochlo...	25021710	1
Fagron	Benzydamine hydrochlo...	13H19-N06	1
Fagron	Benzydamine hydrochlo...	13H19_N06	1
Fagron	Betaine hydrochloride	12C20-N09	1
Fagron	Betaine hydrochloride	24091301	1
Fagron	Betaine hydrochloride	13G29-B02-294870	1
Fagron	Betaine hydrochloride	15121401	1
Fagron	Betaine monohydrate	16H03-B04-330370	1
Fagron	Betamethasone dipropi...	17E22B06336993	1
Fagron	Betamethasone dipropi...	12C21-N13	1
Fagron	Betamethasone dipropi...	17b16b04332218	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Betamethasone dipropi...	18B06B05352731	2
Fagron	Betamethasone dipropi...	18B06-B005-352731	1
Fagron	Betamethasone dipropi...	8K03-N13	1
Fagron	Betamethasone dipropi...	18B06-B05-354807	1
Fagron	Betamethasone dipropi...	17B16-B04-332217	8
Fagron	Betamethasone dipropi...	18b06-b05-350473	1
Fagron	Betamethasone dipropi...	13023004	1
Fagron	Betamethasone dipropi...	18B06B05354807	1
Fagron	Betamethasone dipropi...	17E22-B06-335244	3
Fagron	Betamethasone dipropi...	10E19-N03	1
Fagron	Betamethasone dipropi...	1317-N05	1
Fagron	Betamethasone dipropi...	17B16-B04	1
Fagron	Betamethasone dipropi...	17B16-B04-332224	2
Fagron	Betamethasone dipropi...	17E22-B06-336993	2
Fagron	Betamethasone dipropi...	13I17-N05	2
Fagron	Betamethasone dipropi...	BDP101/0417	1
Fagron	Betamethasone dipropi...	20171127Ma5	1
Fagron	Betamethasone dipropi...	17H16-B04-340537	1
Fagron	Betamethasone dipropi...	17816B04332224	1
Fagron	Betamethasone dipropi...	12B27-N08	3
Fagron	Betamethasone dipropi...	13/17-N05	8
Fagron	Betamethasone dipropi...	BDP014/1115, 2016111	1
Fagron	Betamethasone dipropi...	16C04-B03-320851	1
Fagron	Betamethasone dipropi...	16C04-B03-320852	1
Fagron	Betamethasone dipropi...	16C04-B03-320850	2
Fagron	Betamethasone dipropi...	2109DM0B0011322	1
Fagron	Betamethasone dipropi...	17B16-B04-332219-201	1
Fagron	Betamethasone dipropi...	14C03-B03-296838	1
Fagron	Betamethasone dipropi...	18B06-B05	2
Fagron	Betamethasone dipropi...	11K28-N10	2
Fagron	Betamethasone dipropi...	15H28-B06-313088, 201	1
Fagron	Betamethasone dipropi...	13j16-B01-291844	1
Fagron	Betamethasone dipropi...	9J05-N09	1
Fagron	Betamethasone dipropi...	14C03-B06-296837	7
Fagron	Betamethasone dipropi...	18B06-B05-347338	1
Fagron	Betamethasone dipropi...	15H25-B06-313088	1
Fagron	Betamethasone dipropi...	13J16-B01-291845	1
Fagron	Betamethasone dipropi...	12B27-N08	3
Fagron	Betamethasone dipropi...	12827-N08	1
Fagron	Betamethasone dipropi...	13 17-N050632/1870	1
Fagron	Betamethasone dipropi...	13I17-N05	9
Fagron	Betamethasone dipropi...	18B06B05347338	1
Fagron	Betamethasone dipropi...	11D26-N06	2
Fagron	Betamethasone dipropi...	18B06-B05-346218	4
Fagron	Betamethasone dipropi...	18B06B0534996	1
Fagron	Betamethasone dipropi...	17E22-B06-336705	1
Fagron	Betamethasone dipropi...	13/17-n05	1
Fagron	Betamethasone dipropi...	17B6-B04-332217	1
Fagron	Betamethasone dipropi...	14C03-B06	1
Fagron	Betamethasone dipropi...	16C04-B03-320850, 201	1
Fagron	Betamethasone dipropi...	18B06-B05-350473	3
Fagron	Betamethasone dipropi...	15H28-B06-313088	4
Fagron	Betamethasone valerate	15/16-B07-317428	3
Fagron	Betamethasone valerate	90616-2	1
Fagron	Betamethasone valerate	15F16-B04-315328	3
Fagron	Betamethasone valerate	17329016	2
Fagron	Betamethasone valerate	15/16-B07-317416	3
Fagron	Betamethasone valerate	14H28-B01-305470	1
Fagron	Betamethasone valerate	15F16-B04	2
Fagron	Betamethasone valerate	15F16-804-314054	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Betamethasone valerate	15I16-B07	1
Fagron	Betamethasone valerate	17K-16-B08-347309	2
Fagron	Betamethasone valerate	17G31-B04-344757	2
Fagron	Betamethasone valerate	16B18-B02-320971	7
Fagron	Betamethasone valerate	15I16-B07-319228	2
Fagron	Betamethasone valerate	700103-0004	2
Fagron	Betamethasone valerate	15B24-B04	2
Fagron	Betamethasone valerate	15I16-B07-317416	1
Fagron	Betamethasone valerate	BV106/0717	1
Fagron	Betamethasone valerate	15f16-b04314054	1
Fagron	Betamethasone valerate	17G31B04338575	1
Fagron	Betamethasone valerate	14G31-B04-342639	1
Fagron	Betamethasone valerate	17G31-B04-339784	4
Fagron	Betamethasone valerate	17G31-B04-341576-201	1
Fagron	Betamethasone valerate	17g31-b04-341576	2
Fagron	Betamethasone valerate	17C29B043356575581	1
Fagron	Betamethasone valerate	15824-B04	1
Fagron	Betamethasone valerate	17C29-B04-335654	2
Fagron	Betamethasone valerate	16H22-B06-332351	1
Fagron	Betamethasone valerate	15F16-B04-314054/201	1
Fagron	Betamethasone valerate	17G31-B04-341576	7
Fagron	Betamethasone valerate	17c29-b04-335657	2
Fagron	Betamethasone valerate	15/16-B07-317417	2
Fagron	Betamethasone valerate	17G31-B04	2
Fagron	Betamethasone valerate	17G31-B-04	1
Fagron	Betamethasone valerate	17G31-B04-342639	1
Fagron	Betamethasone valerate	BV-BLT250/0616	1
Fagron	Betamethasone valerate	17K16-B08-345555	2
Fagron	Betamethasone valerate	17G31-B04-339785	1
Fagron	Betamethasone valerate	17K16-B08-345541	7
Fagron	Betamethasone valerate	17G-31-B04-339784	1
Fagron	Betamethasone valerate	17k16-b08-345541	4
Fagron	Betamethasone valerate	18I27-B05	1
Fagron	Betamethasone valerate	14H28-B01-302864/201	1
Fagron	Betamethasone valerate	BV008/0415	1
Fagron	Betamethasone valerate	14H28B01302864	1
Fagron	Betamethasone valerate	17G31-B04-3415765	1
Fagron	Betamethasone valerate	BV009/0117	1
Fagron	Betamethasone valerate	15I16-B07-319832	3
Fagron	Betamethasone valerate	1526260	1
Fagron	Betamethasone valerate	16E23-B03-327354	3
Fagron	Betamethasone valerate	16B18-B02-322109	5
Fagron	Betamethasone valerate	15I16-B07-319832	1
Fagron	Betamethasone valerate	17g31-b04-338575	1
Fagron	Betamethasone valerate	16E23B03328128	1
Fagron	Betamethasone valerate	BV008/0516	1
Fagron	Betamethasone valerate	15/16-B07	1
Fagron	Betamethasone valerate	18G03-B03-351513	1
Fagron	Betamethasone valerate	16H22-B06-330522	1
Fagron	Betamethasone valerate	16H22-B06	2
Fagron	Betamethasone valerate	16H22B06330524	2
Fagron	Betamethasone valerate	18G03-B03-351450	3
Fagron	Betamethasone valerate	17K16-B08-347912	3
Fagron	Betamethasone valerate	17C29-B04-333688	3
Fagron	Betamethasone valerate	17K16-B08	3
Fagron	Betamethasone valerate	16E23-B03-3247990855	1
Fagron	Betamethasone valerate	700103-0002	1
Fagron	Betamethasone valerate	17K16-B08-345542	1
Fagron	Betamethasone valerate	16E23-B03-327350	3
Fagron	Betamethasone valerate	17k16-B08-347309	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Betamethasone valerate	17C29-B04-335676	5
Fagron	Betamethasone valerate	17C29-B04-335657	6
Fagron	Betamethasone valerate	17K16-B08-349755	4
Fagron	Betamethasone valerate	16E23-B03-328128	14
Fagron	Betamethasone valerate	17C29-335654	1
Fagron	Betamethasone valerate	16H22-B06-330523	3
Fagron	Betamethasone valerate	16H22-B06-330524	4
Fagron	Betamethasone valerate	17K16-B08-348538	1
Fagron	Betamethasone valerate	BV008/045	1
Fagron	Betamethasone valerate	17K16-B08-349774	1
Fagron	Betamethasone valerate	16H22-B06-329223	2
Fagron	Betamethasone valerate	17K16-B08-347309	2
Fagron	Betamethasone valerate	16818-B02-320971	1
Fagron	Betamethasone valerate	15F16-B04-314055	1
Fagron	Betamethasone valerate	1516-B07-317416	1
Fagron	Betamethasone valerate	17C29-B04-335658	1
Fagron	Betamethasone valerate	16B18-B02-324798	3
Fagron	Betamethasone valerate	17K16-B08-348042	5
Fagron	Betamethasone valerate	16E23-B03-324800	6
Fagron	Betamethasone valerate	17K16-B08-344754	1
Fagron	Betamethasone valerate	15/16-B07-319832	1
Fagron	Betamethasone valerate	17k16-b08	1
Fagron	Betamethasone valerate	17G31-B04-338575	6
Fagron	Betamethasone valerate	BV036/1215	1
Fagron	Betamethasone valerate	18/27-B05-354880	1
Fagron	Betamethasone valerate	16B18-B02320971	2
Fagron	Betamethasone valerate	16B18-B02-322121	5
Fagron	Betamethasone valerate	17C29-B04-334329	6
Fagron	Betamethasone valerate	B07-317416	1
Fagron	Betamethasone valerate	18I27-B05-356570	1
Fagron	Betamethasone valerate	17C29-B04-336434	4
Fagron	Betamethasone valerate	18G03-B03-351451	1
Fagron	Betamethasone valerate	18/27-B05	1
Fagron	Betamethasone valerate	16818-B02-323033	1
Fagron	Betamethasone valerate	16B18-B02-323033	1
Fagron	Betamethasone valerate	15116-B07	1
Fagron	Betamethasone valerate	16B18-B02-323034	2
Fagron	Betamethasone valerate	16e23-b03-324800	1
Fagron	Betamethasone valerate	18G03-B03-351731	4
Fagron	Betamethasone valerate	16B18B02320971	1
Fagron	Betamethasone valerate	11K-18-N09	1
Fagron	Betamethasone valerate	16B18-B02	4
Fagron	Betamethasone valerate	15116-B07-319832	2
Fagron	Betamethasone valerate	16e23b03324800	1
Fagron	Betamethasone valerate	16E23-B03-324799	10
Fagron	Betamethasone valerate	18G03-B03-352460	2
Fagron	Betamethasone valerate	18I27-B05-354898	1
Fagron	Betamethasone valerate	16b18-b02-324798	1
Fagron	Betamethasone valerate	16H22-B06-329231	4
Fagron	Betamethasone valerate	12L11-N01	1
Fagron	Betamethasone valerate	15F16-B04-314054	5
Fagron	Betamethasone valerate	13K29-B12-288247	1
Fagron	Betamethasone valerate	13K29-B32-28824	1
Fagron	Betamethasone valerate	13K29B32288243	1
Fagron	Betamethasone valerate	14C19-B02-292653	7
Fagron	Betamethasone valerate	13K29-B32-288248	3
Fagron	Betamethasone valerate	14C19-B02-292648	1
Fagron	Betamethasone valerate	13K29-B23-288246	1
Fagron	Betamethasone valerate	NP1351M	1
Fagron	Betamethasone valerate	11K18-N09	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Betamethasone valerate	292652	1
Fagron	Betamethasone valerate	13117-B08-292649	2
Fagron	Betamethasone valerate	13/17-B08-292649	2
Fagron	Betamethasone valerate	14C19-B02-296169	1
Fagron	Betamethasone valerate	14C19-B02	3
Fagron	Betamethasone valerate	14c19-b02-292650	1
Fagron	Betamethasone valerate	12L11-N01	21
Fagron	Betamethasone valerate	14C19-B02-292646	1
Fagron	Betamethasone valerate	13K05-B08-292652	1
Fagron	Betamethasone valerate	13K29-B32-288245	1
Fagron	Betamethasone valerate	13K29-B32-288244	2
Fagron	Betamethasone valerate	14C19-B02-296168	4
Fagron	Betamethasone valerate	13K29-B32-288247	16
Fagron	Betamethasone valerate	13K29-B32-288250	3
Fagron	Betamethasone valerate	2B28-N01	1
Fagron	Betamethasone valerate	13D03-N07	41
Fagron	Betamethasone valerate	13d03-n07	3
Fagron	Betamethasone valerate	NS0706M	3
Fagron	Betamethasone valerate	NS2072M	1
Fagron	Betamethasone valerate	13G19-N01	48
Fagron	Betamethasone valerate	NT0597M	3
Fagron	Betamethasone valerate	13g19-n01	2
Fagron	Betamethasone valerate	13G19N01	1
Fagron	Betamethasone valerate	13G19-N01	1
Fagron	Betamethasone valerate	11E18-N03	2
Fagron	Betamethasone valerate	11A04-N02	1
Fagron	Betamethasone valerate	20101012-1	1
Fagron	Betamethasone valerate	612211181	1
Fagron	Betamethasone valerate	13K29-B32-288246	5
Fagron	Betamethasone valerate	13K29-B32	6
Fagron	Betamethasone valerate	12G09-N07	28
Fagron	Betamethasone valerate	12B28-N01	10
Fagron	Betamethasone valerate	NT059M	1
Fagron	Betamethasone valerate	13K29B32	1
Fagron	Betamethasone valerate	13k29b32288247	1
Fagron	Betamethasone valerate	12K07-N23	8
Fagron	Betamethasone valerate	13K05-B08-292647	2
Fagron	Betamethasone valerate	14H28-B01-302683	1
Fagron	Betamethasone valerate	14G30-B05-309520	6
Fagron	Betamethasone valerate	13K05-B08	1
Fagron	Betamethasone valerate	14H28B01	1
Fagron	Betamethasone valerate	14G30-B05-309519	1
Fagron	Betamethasone valerate	BV-BLT225/0514	4
Fagron	Betamethasone valerate	14G30-B05-309518	2
Fagron	Betamethasone valerate	14G30-B05-29617	1
Fagron	Betamethasone valerate	15B24-B04-307015	3
Fagron	Betamethasone valerate	14c19-b02-296169	1
Fagron	Betamethasone valerate	12K07N23	1
Fagron	Betamethasone valerate	309520	1
Fagron	Betamethasone valerate	14G30-B05-309528	1
Fagron	Betamethasone valerate	15B24-B04-307016	13
Fagron	Betamethasone valerate	14h28-b01-302865	2
Fagron	Betamethasone valerate	14G30-B05-309521	1
Fagron	Betamethasone valerate	14H28-B01-302864	6
Fagron	Betamethasone valerate	14H28-B01-302877	4
Fagron	Betamethasone valerate	14H28-B01-302866	4
Fagron	Betamethasone valerate	14G30-B05	4
Fagron	Betamethasone valerate	14G30B05309519	1
Fagron	Betamethasone valerate	14H28-B01-302865	13
Fagron	Betamethasone valerate	14G30B05296172	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Betamethasone valerate	14G30-B05-296172	3
Fagron	Betamethasone valerate	14G30-B05-296173	11
Fagron	Betamethasone valerate	14H28-B01	3
Fagron	Betamethasone valerate	14G30-Bo5-296172	1
Fagron	Betamethasone valerate	14C19B02	1
Fagron	Betamethasone valerate	14G30-B05-296175	2
Fagron	Betamethasone valerate	13K05-B08-291663	1
Fagron	Betamethasone valerate	BV019/0514	3
Fagron	Betamethasone, micron...	14C13-B02-296154	1
Fagron	Betamethasone, micron...	14C13-B02	2
Fagron	Betamethasone, micron...	14C13-B02-305561	1
Fagron	Betamethasone, micron...	14C13-B02-296164	1
Fagron	Betamethasone, micron...	15L09-B05_322119	1
Fagron	Betamethasone, micron...	9K25-N13	1
Fagron	Betamethasone, micron...	9J26-N08	1
Fagron	Betamethasone, micron...	12J01-N02	1
Fagron	Betamethasone, micron...	12D26-N08	1
Fagron	Betamethasone, micron...	11H31-B06-288217	1
Fagron	Betamethasone, micron...	14c13-bo2-305561	1
Fagron	Betamethasone, micron...	14C13-B02-296164	1
Fagron	Betamethasone, micron...	9071304	1
Fagron	Betamethasone, micron...	15L09-B05	1
Fagron	Betamethasone, micron...	16I07-B01-341871	2
Fagron	Betamethasone, micron...	12D26-N08	1
Fagron	Betamethasone-17-vale...	12G02-N02	2
Fagron	Betamethasone-17-vale...	12L19-N06	4
Fagron	Betamethasone-17-vale...	12C19-N07	1
Fagron	Betamethasone-17-vale...	12J19-N06	1
Fagron	Betamethasone-17-vale...	12G02-N15	3
Fagron	Betamethasone-17-vale...	13E08-N01	1
Fagron	Bifonazole	16J03-B08-331029	2
Fagron	Bifonazole	18A15-B05_344514	1
Fagron	Bifonazole	16J03-B08-333936	2
Fagron	Bifonazole	16D07-B10-321576	1
Fagron	Bifonazole	16D07-B10-324381	1
Fagron	Bifonazole	16J03-B08-327046	2
Fagron	Bifonazole	18A15-B05-344515	1
Fagron	Bifonazole	12E29-N03	3
Fagron	Bifonazole	1607-B10-324380	1
Fagron	Bifonazole	16D07-B10-32480	1
Fagron	Bifonazole	14K10-B01-316618	1
Fagron	Bifonazole	12E29-N03	2
Fagron	Bifonazole	14K10-B01-312067	2
Fagron	Bifonazole	14K10-B01-308349	2
Fagron	Bifonazole	14K10-B01-306783	2
Fagron	Bifonazole	14K10-B01-300541	6
Fagron	Bifonazole	14C28-B03-292388	1
Fagron	Bifonazole	13G10-N02	1
Fagron	Bifonazole	13G10-N02	1
Fagron	Bifonazole	18A15-B05-350644	1
Fagron	Bifonazole	11F20-N06	1
Fagron	Bifonazole	16D07-B10-324380	1
Fagron	Bifonazole	18A15-B05-357781	1
Fagron	Bifonazole	14K10-B01-306783	2
Fagron	Bifonazole	12G11-N09	2
Fagron	Biotin	17K20-B03-351524	1
Fagron	Biotin	14d04-b07-292963	1
Fagron	Biotin	13K29-N02	2
Fagron	Biotin	20110634	1
Fagron	Biotin	12111302	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Biotin	12C14-N01	1
Fagron	Biotin	13081203	1
Fagron	Biotin	13111202	1
Fagron	Birch leaf, cut	7061301	1
Fagron	Bismuth gallath, basic	15J26-B02-320473	6
Fagron	Bismuth gallath, basic	15J26-B02-317883	1
Fagron	Bismuth gallath, basic	1318B1329409	1
Fagron	Bismuth gallath, basic	15J26-B02-317883	6
Fagron	Bismuth gallath, basic	2015000361	2
Fagron	Bismuth gallath, basic	15H05-B07	1
Fagron	Bismuth gallath, basic	2013006220	1
Fagron	Bismuth gallath, basic	2015000866	1
Fagron	Bismuth gallath, basic	15C27-B02-311879	1
Fagron	Bismuth gallath, basic	15b02-B04-309481	1
Fagron	Bismuth gallath, basic	15B02-B04-311878	1
Fagron	Bismuth gallath, basic	15B02-B04	2
Fagron	Bismuth gallath, basic	15B02-B04-304408	6
Fagron	Bismuth gallath, basic	14H13-B07-303639	5
Fagron	Bismuth gallath, basic	14H13-B07	1
Fagron	Bismuth gallath, basic	15B02-B04-304413	5
Fagron	Bismuth gallath, basic	15J26-B02-321313	7
Fagron	Bismuth gallath, basic	12B02-B04-304408	1
Fagron	Bismuth gallath, basic	14H13B07303639	1
Fagron	Bismuth gallath, basic	14H13-B07-303639	1
Fagron	Bismuth gallath, basic	15B02-B04-304408	1
Fagron	Bismuth gallath, basic	13I18-B13-294108	6
Fagron	Bismuth gallath, basic	13i18-b13-294108	1
Fagron	Bismuth gallath, basic	1318B13294108	1
Fagron	Bismuth gallath, basic	13I18-B13-294098	2
Fagron	Bismuth gallath, basic	2014006082	1
Fagron	Bismuth gallath, basic	15B02-B04-209481	1
Fagron	Bismuth gallath, basic	15H05-B07-311880	5
Fagron	Bismuth gallath, basic	2015007568	4
Fagron	Bismuth gallath, basic	1726B01341866	1
Fagron	Bismuth gallath, basic	17I26-B01-341866	5
Fagron	Bismuth gallath, basic	17I26-B01-340552	1
Fagron	Bismuth gallath, basic	13/18-B13-292057	1
Fagron	Bismuth gallath, basic	18A29-B06-346303	1
Fagron	Bismuth gallath, basic	18A29-B06-3460303	1
Fagron	Bismuth gallath, basic	180808186	1
Fagron	Bismuth gallath, basic	2018000237	1
Fagron	Bismuth gallath, basic	18A29B06-346305	1
Fagron	Bismuth gallath, basic	18A29-B06-346305	1
Fagron	Bismuth gallath, basic	18G06-B10-352463	1
Fagron	Bismuth gallath, basic	17G06-B10-355772	1
Fagron	Bismuth gallath, basic	18G06-B10-353902	1
Fagron	Bismuth gallath, basic	19B05-B09-191022	1
Fagron	Bismuth gallath, basic	18G06-B10-353916	3
Fagron	Bismuth gallath, basic	16K17-B04-3355593	1
Fagron	Bismuth gallath, basic	2017006945	3
Fagron	Bismuth gallath, basic	180307BB/2016006984	1
Fagron	Bismuth gallath, basic	15J26-B02	3
Fagron	Bismuth gallath, basic	15/26-B02-320473	1
Fagron	Bismuth gallath, basic	15126-B02-321313	1
Fagron	Bismuth gallath, basic	15J26-B02-325472	5
Fagron	Bismuth gallath, basic	325472	1
Fagron	Bismuth gallath, basic	15J26-B022	1
Fagron	Bismuth gallath, basic	15J26-B02-321313	2
Fagron	Bismuth gallath, basic	16K17-B04-335597	11
Fagron	Bismuth gallath, basic	16K17B04335597	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Bismuth gallath, basic	16K17B04	1
Fagron	Bismuth gallath, basic	16K17-B04	1
Fagron	Bismuth gallath, basic	17E16-B02	3
Fagron	Bismuth gallath, basic	17/26-B01-340556	4
Fagron	Bismuth gallath, basic	2016006984	1
Fagron	Bismuth gallath, basic	17I26-B01-340556	3
Fagron	Bismuth gallath, basic	13I18-b13-292057	1
Fagron	Bismuth gallath, basic	13/18-B13294098	1
Fagron	Bismuth gallath, basic	15b02b04304408	1
Fagron	Bismuth gallath, basic	13/18-B-13-292057	1
Fagron	Bismuth gallath, basic	12E16-N03	1
Fagron	Bismuth gallath, basic	2012006724	2
Fagron	Bismuth gallath, basic	12K26-N01	2
Fagron	Bismuth gallath, basic	11F15-N03	1
Fagron	Bismuth gallath, basic	13E21-N03	1
Fagron	Bismuth gallath, basic	2010002685	2
Fagron	Bismuth gallath, basic	411092-3-P	1
Fagron	Bismuth gallath, basic	13F10-N12	1
Fagron	Bismuth gallath, basic	13G10-N12	1
Fagron	Bismuth gallath, basic	10H31-N01	1
Fagron	Bismuth gallath, basic	13K08-N03	2
Fagron	Bismuth gallath, basic	13/18-B13-294098	2
Fagron	Bismuth gallath, basic	13G10-N12	3
Fagron	Bismuth gallath, basic	13D19-N02	10
Fagron	Bismuth gallath, basic	2012006724	2
Fagron	Bismuth gallath, basic	13/18-B13-292058	3
Fagron	Bismuth gallath, basic	1318-B13-292057	2
Fagron	Bismuth gallath, basic	13/18-B13-294108	3
Fagron	Bismuth gallath, basic	13I18-B13-294098	2
Fagron	Bismuth gallath, basic	13 18-B13-294098	1
Fagron	Bismuth gallath, basic	13179111	3
Fagron	Bismuth gallath, basic	13I18-B13	3
Fagron	Bismuth gallath, basic	13G10-N12	1
Fagron	Bismuth gallath, basic	13I18-B13	3
Fagron	Bismuth gallath, basic	12K26-N01	1
Fagron	Bismuth gallath, basic	1318b13292058	1
Fagron	Bismuth gallath, basic	13K08-N03	1
Fagron	Bismuth gallath, basic	5121	1
Fagron	Bismuth gallath, basic	13I18-B13-292058	1
Fagron	Bismuth gallath, basic	13I18-B13-292057	1
Fagron	Bismuth gallath, basic	13K08-N03	1
Fagron	Bismuth gallath, basic	12/13-N11	2
Fagron	Bitter fennel		1
Fagron	Bitter fennel	1406211545	1
Fagron	Bitter fennel	6858913	1
Fagron	Boric acid	11d13-n08	1
Fagron	Boric acid	12G04-N10	1
Fagron	Boric acid	14133101	1
Fagron	Boric acid	13I04-N06	2
Fagron	Boric acid	11I26-n05	1
Fagron	Boric acid	14B07-B02-295752	1
Fagron	Boric acid	12f01-n08	1
Fagron	Boric acid	9082012B	1
Fagron	Budesonide, micronized	13H22-N03	1
Fagron	Budesonide, micronized	18172601	1
Fagron	Budesonide, micronized	12G23-N01	2
Fagron	Budesonide, micronized	15E27_B10	1
Fagron	Budesonide, micronized	15E27-B10-308322	1
Fagron	Budesonide, micronized	17A17-B01-331248	1
Fagron	Budesonide, micronized	14H15-B02-301972	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Budesonide, micronized	15E27-B01-308321	1
Fagron	Budesonide, micronized	14H15-B02-301971	2
Fagron	Budesonide, micronized	14H15-B02-297926	1
Fagron	Budesonide, micronized	17A17-B01-331448	1
Fagron	Budesonide, micronized	14A30-B09_291366	1
Fagron	Budesonide, micronized	13A07-N09	1
Fagron	Budesonide, micronized	14H15-B02	3
Fagron	Budesonide, micronized	14A30-B09-291366	2
Fagron	Budesonide, micronized	15E27-B10-315704	1
Fagron	Budesonide, micronized	17F28-B04-345016	1
Fagron	Butylated hydroxytolu...	11308701	1
Fagron	Butylated hydroxytolu...	15426801	11
Fagron	Caffeine	15H11-B02	1
Fagron	Caffeine citrate	2813E1282	1
Fagron	Caffeine citrate	11082015A	1
Fagron	Caffeine citrate	15A05-B03-307135	3
Fagron	Caffeine citrate	3914U-01282	1
Fagron	Caffeine citrate	15A05-B03-318006	13
Fagron	Caffeine citrate	15a05-b03-318006	1
Fagron	Caffeine citrate	15a05-b03-302746	1
Fagron	Caffeine citrate	13B08N02	1
Fagron	Caffeine citrate	15A05B03	1
Fagron	Caffeine citrate	15A05-B03	3
Fagron	Caffeine citrate	15A05-B03	1
Fagron	Caffeine citrate	3513I-01282	1
Fagron	Caffeine citrate	13D30-N02	1
Fagron	Caffeine citrate	11C29-N04	1
Fagron	Caffeine citrate	12C16-N03	1
Fagron	Caffeine citrate	204311282	1
Fagron	Caffeine citrate	4213U01282	1
Fagron	Caffeine citrate	11C29-N04	1
Fagron	Caffeine citrate	12J05-N06	1
Fagron	Caffeine citrate	10H02-N02	2
Fagron	Caffeine citrate	5013U01282	1
Fagron	Caffeine citrate	2002311282	1
Fagron	Caffeine citrate	26021401	1
Fagron	Caffeine citrate	CCC1205	1
Fagron	Caffeine citrate	3314E-01282	1
Fagron	Caffeine citrate	2714E-01282	1
Fagron	Caffeine citrate	13D30-No2	1
Fagron	Caffeine citrate	CCD-1401	3
Fagron	Calcium citrate tetra...	13d16-n06	1
Fagron	Calcium citrate tetra...	13D16-N06	1
Fagron	Calcium citrate tetra...	150903901	1
Fagron	Calcium gluconate	12D17-N01	1
Fagron	Calcium gluconate	15D13-B98-307445	1
Fagron	Calcium gluconate	14D11-B04	1
Fagron	Calcium gluconate	14D11-B04-294136	2
Fagron	Calcium gluconate	10G22-N12	1
Fagron	Calcium glycerophosph...	1317-N04	2
Fagron	Calcium glycerophosph...	1216000005	1
Fagron	Calcium glycerophosph...	11G19-N05	1
Fagron	Calcium glycerophosph...	12I18-N06	3
Fagron	Calcium glycerophosph...	13A18-N02	5
Fagron	Calcium glycerophosph...	3115A-01212	1
Fagron	Calcium glycerophosph...	14B10-B02	1
Fagron	Calcium glycerophosph...	15B03B03	1
Fagron	Calcium glycerophosph...	5212M-01212	1
Fagron	Calcium glycerophosph...	14B10B02	1
Fagron	Calcium glycerophosph...	14B10-B02-296841	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Calcium glycerophosph...	14B10-B02-291642	1
Fagron	Calcium glycerophosph...	1314000019	1
Fagron	Calcium glycerophosph...	15B03-B03	1
Fagron	Calcium glycerophosph...	13A18N02	2
Fagron	Calcium lactate penta...	14H01-B05-334588	1
Fagron	Calcium lactate penta...	12/12-N06	1
Fagron	Calcium lactate penta...	12I12-N06	2
Fagron	Calcium lactate penta...	13L18-B04-295541	2
Fagron	Calcium lactate penta...	14H01-B05	1
Fagron	Calcium lactate penta...	13L18B04	1
Fagron	Calcium lactate penta...	17J18-B01-347174	1
Fagron	Calcium lactate penta...	13L18-B04-295541	1
Fagron	Calcium lactate penta...	12I12-N01	1
Fagron	Calcium pantothenate	14G02-B10-297373	1
Fagron	Calcium pantothenate	13H28-N11	1
Fagron	Calcium pantothenate	12F19-N07	1
Fagron	Calcium pantothenate	1111110405	1
Fagron	Calendula flower, cut	1241B120312	4
Fagron	Camphor	18066404	4
Fagron	Camphor	27011	1
Fagron	Camphor	15/18-B11-321658	1
Fagron	Camphor	13K28-B04-289340	1
Fagron	Camphor	1309019-02	1
Fagron	Camphor	15097801	1
Fagron	Camphor	14G23-B04-301879	1
Fagron	Camphor	14G23-B04-301878	2
Fagron	Camphor	13K28-B04-294937	1
Fagron	Camphor	14G23-B04	2
Fagron	Camphor	13K28_B04_294933	1
Fagron	Camphor	13K28-B04	1
Fagron	Camphor	12B14-N09	1
Fagron	Camphor	11H30-N15	1
Fagron	Camphor	12D23-N03	1
Fagron	Camphor	11K23-N04	1
Fagron	Camphor	690442	1
Fagron	Camphor	14G23-B04-299202	2
Fagron	Camphor	13H23-B04-288376	1
Fagron	Camphor	13H23-B40-291888	1
Fagron	Camphor	12G02-N09	1
Fagron	Camphor	13H23-B40-294116	1
Fagron	Camphor	8082014A	1
Fagron	Camphor	13h23-b40-294116	1
Fagron	Camphor	825653	1
Fagron	Camphor	13H23-B40-296252	1
Fagron	Camphor	13K28-B04	2
Fagron	Camphor	14G23-B04-296868	1
Fagron	Capsaicin	CP-12003	1
Fagron	Captopril	17G27-B02-337989	3
Fagron	Captopril	80616-6	1
Fagron	Captopril	14B26-B40-318894	6
Fagron	Captopril	14B26-B01-314460	1
Fagron	Captopril	700766	1
Fagron	Captopril	14B26-B40-320183	8
Fagron	Captopril	14B26-B40-321319	5
Fagron	Captopril	17g27-b02-337989	1
Fagron	Captopril	14B26-B01-302337	1
Fagron	Captopril	14B26-B40-325275	5
Fagron	Captopril	14B26-B01-312221	1
Fagron	Captopril	EC140101	1
Fagron	Captopril	14B26-B01-291555	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Captopril	14B26B01	1
Fagron	Captopril	14B26-B01_302338	1
Fagron	Captopril	14B26-B01-291555	2
Fagron	Captopril	14B26-B01-291556	2
Fagron	Captopril	14B26-W0129-1557	1
Fagron	Captopril	12G25-N08	4
Fagron	Captopril	14B26-B01_2911557	1
Fagron	Captopril	12K-16N04	1
Fagron	Captopril	14B26-B01	4
Fagron	Captopril	14B26-B01-291557	5
Fagron	Captopril	17G27-B02-340798	2
Fagron	Captopril	12K16-N04	1
Fagron	Captopril	14B26-B01-302338	2
Fagron	Captopril	17G27-B02-346203	2
Fagron	Captopril	14B26-B01-316257	2
Fagron	Captopril	14196902	1
Fagron	Captopril	17G27-B02-352730	1
Fagron	Captopril	17G27-B02-356147	1
Fagron	Captopril	17G-B02-352730	1
Fagron	Carmellose sodium (40...	13D22-N02	1
Fagron	Carmellose sodium (40...	15G03-Bo2-312862	1
Fagron	Cellulose, microcryst...	17A11-B04-331463	1
Fagron	Cellulose, microcryst...	15C03-B06-305242	1
Fagron	Cellulose, microcryst...	2013M-01350	1
Fagron	Cellulose, microcryst...	9K10-N10	1
Fagron	Cellulose, microcryst...	13A03-N01	1
Fagron	Cellulose, microcryst...	17A11-B04-340541	1
Fagron	Cetostearyl alcohol	13L17-B05	1
Fagron	Cetostearyl alcohol	16B23-B07-318218	2
Fagron	Cetostearyl alcohol	13L17-B05-289805	1
Fagron	Cetostearyl alcohol	8256792	2
Fagron	Cetostearyl alcohol	13F07-N05	2
Fagron	Cetostearyl alcohol	15D24-B02-311993	2
Fagron	Cetostearyl alcohol	15C31-B99-307062	2
Fagron	Cetostearyl alcohol	13L17-B05-289806	2
Fagron	Cetostearyl alcohol	15D24-B02-308261	1
Fagron	Cetostearyl alcohol	14G31-B03-297611	1
Fagron	Cetostearyl alcohol	15D24-B02	1
Fagron	Cetostearyl alcohol	12C30-N01	1
Fagron	Cetostearyl alcohol	9652865	1
Fagron	Cetostearyl alcohol	13C13-N13	1
Fagron	Cetostearyl alcohol	12J29-N04	1
Fagron	Cetostearyl alcohol	12L19-N03	2
Fagron	Cetostearyl alcohol	11E17-N13	1
Fagron	Cetostearyl alcohol	12K22-N03	3
Fagron	Cetostearyl alcohol	12/17-N05	1
Fagron	Cetostearyl alcohol	12J17-N05	2
Fagron	Cetostearyl alcohol	12E30-N13	1
Fagron	Cetostearyl alcohol	3414E-01078	1
Fagron	Cetostearyl alcohol	16B23-B07-320384	1
Fagron	Cetostearyl alcohol	12J17-N05	1
Fagron	Cetostearyl alcohol	17F21-B01-341021	1
Fagron	Cetostearyl alcohol	13L17-B05-289806	1
Fagron	Cetostearyl alcohol	7886542	1
Fagron	Cetostearyl alcohol	12G02-N04	7
Fagron	Cetostearyl alcohol	12g02-n04	1
Fagron	Cetostearyl alcohol	11f30-n11	1
Fagron	Cetostearyl alcohol	10F17-N02	1
Fagron	Cetostearyl alcohol	17F21-B01	1
Fagron	Cetostearyl alcohol	13L17-B05-294999	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Cetyl alcohol	15D24-B02-311993	1
Fagron	Cetyl alcohol	12A-02N12	1
Fagron	Cetyl alcohol	18A23-B04-349155	1
Fagron	Cetyl palmitate 15	12C27-N07	2
Fagron	Cetyl palmitate 15	12D19-N02	1
Fagron	Chloral hydrate	160609/1	1
Fagron	Chloral hydrate	12H06N08	2
Fagron	Chloral hydrate	12C26-M04	1
Fagron	Chloral hydrate	11J05-N01	1
Fagron	Chloral hydrate	11D13-N07	1
Fagron	Chloral hydrate	12H06-N080101/1612	1
Fagron	Chloral hydrate	13F06-N070631/786	1
Fagron	Chloral hydrate	12H06-N08	9
Fagron	Chloral hydrate	13K05-B13	1
Fagron	Chloral hydrate	12H06-N080160/161	1
Fagron	Chloral hydrate	12C26-M05	1
Fagron	Chloral hydrate	12H06N8	1
Fagron	Chloral hydrate	14218710	2
Fagron	Chloral hydrate	12H06_N08	1
Fagron	Chloral hydrate	13D22-N11	1
Fagron	Chloramphenicol	16/12-B07	1
Fagron	Chloramphenicol	11J14-N03	1
Fagron	Chloramphenicol	12L11-N06	9
Fagron	Chloramphenicol	12B13-N06	6
Fagron	Chloramphenicol	13A11N02	1
Fagron	Chloramphenicol	11/14-N03	2
Fagron	Chloramphenicol	17A23-B-08-336728	1
Fagron	Chloramphenicol	13A11-N02	9
Fagron	Chloramphenicol	12F01-N23	1
Fagron	Chloramphenicol	14E01-B02-296625	1
Fagron	Chloramphenicol	15G14-B04-322516	4
Fagron	Chloramphenicol	Cealo/Phönix	1
Fagron	Chloramphenicol	14E01-B02-308559	5
Fagron	Chloramphenicol	14E01-B02-306146	5
Fagron	Chloramphenicol	14E01-B02-296623	2
Fagron	Chloramphenicol	14E01-B02	4
Fagron	Chloramphenicol	DY00912129	1
Fagron	Chloramphenicol	2111I-01259	1
Fagron	Chloramphenicol	13A31-B02-290263	1
Fagron	Chloramphenicol	16/12-B07-332100	1
Fagron	Chloramphenicol	13A31-B02-290300	3
Fagron	Chloramphenicol	13A31-B02_286708	1
Fagron	Chloramphenicol	13A31-B02-286709	3
Fagron	Chloramphenicol	16I12-B07-326037	1
Fagron	Chloramphenicol	17A23-B08-336728	1
Fagron	Chloramphenicol	13A31-B02-286708	4
Fagron	Chloramphenicol	14E01-B02-296624	3
Fagron	Chloramphenicol	14E01B02296624	3
Fagron	Chloramphenicol	1202500003	1
Fagron	Chloramphenicol	14e01-b02-296624	1
Fagron	Chloramphenicol	Ch.B-A250712	1
Fagron	Chloramphenicol	12A11-N02	1
Fagron	Chloramphenicol	17A23B08337138	1
Fagron	Chloramphenicol	16112-B07-326036	1
Fagron	Chlorhexidine diacetate	12E14-N11	1
Fagron	Chlorhexidine diacetate	15J02-B08	1
Fagron	Chlorhexidine diacetate	A-130037	1
Fagron	Chlorhexidine diacetate	1-8991-1-01-01	1
Fagron	Chlorhexidine diacetate	15J02-B08-322916	2
Fagron	Chlorhexidine diacetate	13E03-B06-287926	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Chlorhexidine diacetate	13E03-B06	1
Fagron	Chlorhexidine diacetate	15438901	1
Fagron	Chlorhexidine dihydro...	14A09-B03-293425	1
Fagron	Chlortetracycline hyd...	14J08-B07-300250	5
Fagron	Chlortetracycline hyd...	14J08-B07-303350	1
Fagron	Chlortetracycline hyd...	16B10-B09-318088	5
Fagron	Chlortetracycline hyd...	16B10-B09-322715	2
Fagron	Chlortetracycline hyd...	FC11051038	1
Fagron	Chlortetracycline hyd...	12K30-N07	1
Fagron	Chlortetracycline hyd...	13111-N07	1
Fagron	Chlortetracycline hyd...	13/11-N07	2
Fagron	Chlortetracycline hyd...	14J08-B07-307100	4
Fagron	Chlortetracycline hyd...	16B10-B09-326865	1
Fagron	Chlortetracycline hyd...	16B10-B09-331443	1
Fagron	Chlortetracycline hyd...	13/11-N07	2
Fagron	Chlortetracycline hyd...	16B10-B09-336332	3
Fagron	Chlortetracycline hyd...	16b10-b09-336332	1
Fagron	Chlortetracycline hyd...	16b10-b09	1
Fagron	Chlortetracycline hyd...	14J08-B07-303352	5
Fagron	Chlortetracycline hyd...	16B10-B09-342410	1
Fagron	Chlortetracycline hyd...	FC13030134	1
Fagron	Chlortetracycline hyd...	16B10-B09-351565	1
Fagron	Chlortetracycline hyd...	16B10-B09-352970	1
Fagron	Chlortetracycline hyd...	16B10-B09-348359	1
Fagron	Chlortetracycline hyd...	27121306	1
Fagron	Chlortetracycline hyd...	12H28-N09	4
Fagron	Chlortetracycline hyd...	12K30-N07	9
Fagron	Chlortetracycline hyd...	12H28N09	1
Fagron	Chlortetracycline hyd...	12k30-N07	1
Fagron	Chlortetracycline hyd...	14J08-B07-300250	2
Fagron	Chlortetracycline hyd...	14J08-B07	4
Fagron	Chlortetracycline hyd...	13K26-N06	1
Fagron	Chlortetracycline hyd...	14J08-B07-307100	2
Fagron	Chlortetracycline hyd...	14J08B07307100	3
Fagron	Cholesterol	15G31-B05-314182	2
Fagron	Cholesterol	31.05.2018	1
Fagron	Cholesterol	15G31-B05-316721	2
Fagron	Cholesterol	15B09	1
Fagron	Cholesterol	15J23-B03-320483	5
Fagron	Cholesterol	15J23-B03-325106	3
Fagron	Cholesterol	15J23-B03-328249	2
Fagron	Cholesterol	15J23-B03-322902	1
Fagron	Cholesterol	15J23B03320483	1
Fagron	Cholesterol	16G05-B03-331396	2
Fagron	Cholesterol	17E17-B05-339074	4
Fagron	Cholesterol	15G31-B05_316721	1
Fagron	Cholesterol	17K16-B11-343922	2
Fagron	Cholesterol	16G05-B03-334366	1
Fagron	Cholesterol	17K16-B11-348761	1
Fagron	Cholesterol	17E17-B05-337588	2
Fagron	Cholesterol	17K16-B11-350557	1
Fagron	Cholesterol	15G31-B05-312377	1
Fagron	Cholesterol	15B09-B02-308259	4
Fagron	Cholesterol	8J28-N10	1
Fagron	Cholesterol	11A17-N07	1
Fagron	Cholesterol	13F10-B40-289016	2
Fagron	Cholesterol	10B15-N07	1
Fagron	Cholesterol	14B27-B05-293548	2
Fagron	Cholesterol	13F10-B40	1
Fagron	Cholesterol	14B27-B05	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Cholesterol	14B27-B05-293550	1
Fagron	Cholesterol	14G04-B02-300830	2
Fagron	Cholesterol	13G10-N11	3
Fagron	Cholesterol	14G04-B02-300831	1
Fagron	Cholesterol	15G31-B05-312378	1
Fagron	Cholesterol	15B09-B02	1
Fagron	Cholesterol	12I13-N02	2
Fagron	Cholesterol	13C25-N12	2
Fagron	Cholesterol	2013U-01270	1
Fagron	Cholesterol	10H24-N06	1
Fagron	Cholesterol	15G31-B05	2
Fagron	Cholesterol	14G04-B02	1
Fagron	Cholesterol	12D02-N10	1
Fagron	Chondroitin sulfate ...	1821031801	1
Fagron	Ciclopirox olamine	10L09-N03	1
Fagron	Ciclopirox olamine	14F13-B03-302339	1
Fagron	Ciclopirox olamine	14E21-M05	1
Fagron	Ciclopirox olamine	11h03-n02	1
Fagron	Ciclopirox olamine	13A16-N04	2
Fagron	Ciclopirox olamine	12G12-N09	2
Fagron	Ciclopirox olamine	A160013	1
Fagron	Ciclosporin	14a15-n01	1
Fagron	Ciclosporin	15f29-b14-316270	2
Fagron	Ciclosporin	13H19-N09	1
Fagron	Ciclosporin	15D23-B05-307217	1
Fagron	Ciclosporin	14718-b02-300180	1
Fagron	Ciclosporin	15D23-B05	2
Fagron	Ciclosporin	13B19N02	1
Fagron	Ciclosporin	13A28-N04	1
Fagron	Ciclosporin	13H19-N09	2
Fagron	Ciclosporin	14/18-B02-300180	1
Fagron	Ciclosporin	14A15-N01	3
Fagron	Ciclosporin	14G22-B02-296891	5
Fagron	Ciclosporin	15F29B14316270	1
Fagron	Ciclosporin	13B19-N02	2
Fagron	Ciclosporin	13H19-N09	2
Fagron	Ciclosporin	13414038	1
Fagron	Ciclosporin	15D23-B05-308853	5
Fagron	Ciclosporin	1103531	1
Fagron	Ciclosporin	15F29-B14-316270	1
Fagron	Ciclosporin	13B19-N02	1
Fagron	Ciclosporin	12E31-N02	7
Fagron	Ciclosporin	11K17-N13	1
Fagron	Ciclosporin	13D26-N07	5
Fagron	Ciclosporin	13A28-n04	1
Fagron	Ciclosporin	13A28-N04	6
Fagron	Ciclosporin	13D26-N07	5
Fagron	Ciclosporin	13A28-N04	1
Fagron	Ciclosporin	13B19N02	1
Fagron	Ciclosporin	10J18-N06	1
Fagron	Ciclosporin	13417083	1
Fagron	Ciclosporin	18G30-B04-354759	1
Fagron	Ciclosporin	18A18-B07-344620	1
Fagron	Ciclosporin	18136104	4
Fagron	Ciclosporin	15F29-B14	1
Fagron	Ciclosporin	15f29-b14-325654	1
Fagron	Ciclosporin	25051302	1
Fagron	Cinchocaine hydrochlo...	12F25-N10	1
Fagron	Cinnarizine	8D11-N04	1
Fagron	Cinnarizine	14D09-B01-292883	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Cinnarizine	13H01-N07	4
Fagron	Cinnarizine	13I26-N22	1
Fagron	Cinnarizine	14D09-B01	3
Fagron	Cinnarizine	13I26-N22	1
Fagron	Ciprofloxacin hydroch...	14I15-B08	1
Fagron	Ciprofloxacin hydroch...	14I15-B08	1
Fagron	Ciprofloxacin hydroch...	15D07-B01-312864	1
Fagron	Ciprofloxacin hydroch...	14B25-B02-291-783	1
Fagron	Ciprofloxacin hydroch...	14I15-B08-303363	1
Fagron	Ciprofloxacin hydroch...	14/15-B08-303363	1
Fagron	Ciprofloxacin hydroch...	14B25-B02-291783	2
Fagron	Citric acid	14G04B09293876	2
Fagron	Citric acid	14K05_B04_302815	1
Fagron	Citric acid	15g16-b01-330450	2
Fagron	Citric acid	330449	1
Fagron	Citric acid	15G16-B01-330450	10
Fagron	Citric acid	13G24-N240739/900	1
Fagron	Citric acid	28021404	1
Fagron	Citric acid	14G04-b09-293876	1
Fagron	Citric acid	15G16-B01330450	1
Fagron	Citric acid	14G04-B09-293874	12
Fagron	Citric acid	13G24-N24	12
Fagron	Citric acid	14K05-B04-302817	2
Fagron	Citric acid	1181765	4
Fagron	Citric acid	15G16B01330450	1
Fagron	Citric acid	UJ501	2
Fagron	Citric acid	15I01-B06-324348	1
Fagron	Citric acid	11K11-N020444/1975	1
Fagron	Citric acid	13K18-B34-288373	1
Fagron	Citric acid	11C24N022869/4250	1
Fagron	Citric acid	13k18-b34-291522	1
Fagron	Citric acid	14E06-B07-293875	3
Fagron	Citric acid	13H22-N09	2
Fagron	Citric acid	13B04-N01	1
Fagron	Citric acid	14E06-B07293872	1
Fagron	Citric acid	18F25-B02-355298	1
Fagron	Citric acid	351453	1
Fagron	Citric acid	17E15-B04-347169	3
Fagron	Citric acid	Ch.B. : 18025702	1
Fagron	Citric acid	17E15-B04-345678	1
Fagron	Citric acid	17E15-B04	2
Fagron	Citric acid	15G16b01	1
Fagron	Citric acid	15G16-B01330449	2
Fagron	Citric acid	15G16B01	1
Fagron	Citric acid	15G16-B01	2
Fagron	Citric acid	14b06-b07-293872	1
Fagron	Citric acid	10G21-N05	2
Fagron	Citric acid	14G04-B09-309836	3
Fagron	Citric acid	12B01-N01	1
Fagron	Citric acid	13A04-N04	1
Fagron	Citric acid	13A-04-N04	2
Fagron	Citric acid	12B01-No1	1
Fagron	Citric acid	13A04-N040019/4900	1
Fagron	Citric acid	14E-06-B07-293890	1
Fagron	Citric acid	13A04-N04	20
Fagron	Citric acid	14k05-b04-302816	1
Fagron	Citric acid	12B01N01	1
Fagron	Citric acid	13A04N04	1
Fagron	Citric acid	14G04-B09-309707	1
Fagron	Citric acid	13G04-N05	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Citric acid	26082013B	1
Fagron	Citric acid	12801-N01	1
Fagron	Citric acid	12b01-N01	2
Fagron	Citric acid	14E06-B07-293872	4
Fagron	Citric acid	13804-N02	1
Fagron	Citric acid	12F06-M06	3
Fagron	Citric acid	13K18-B34_291522	1
Fagron	Citric acid	11C24-N02	5
Fagron	Citric acid	11K07-N01	10
Fagron	Citric acid	20071203	1
Fagron	Citric acid	700309-0003	2
Fagron	Citric acid	13L05B02296247	1
Fagron	Citric acid	1913A-01024	1
Fagron	Citric acid	700309-001	1
Fagron	Citric acid	1151608	2
Fagron	Citric acid	12B01-N01	31
Fagron	Citric acid	4810E-01024	1
Fagron	Citric acid	12B01-N01	1
Fagron	Citric acid	1155609	13
Fagron	Citric acid	probe	1
Fagron	Citric acid	11K07N01	1
Fagron	Citric acid	13k18-B34-288372	1
Fagron	Citric acid	14G04-B09-293877	9
Fagron	Citric acid	14E06-B07-293890	3
Fagron	Citric acid	16072014B	1
Fagron	Citric acid	13K18-B34-291072	3
Fagron	Citric acid	13F12-N01	3
Fagron	Citric acid	13L05-B02-309050	1
Fagron	Citric acid	13L05-B02-296246	3
Fagron	Citric acid	13K18-B34-291521	4
Fagron	Citric acid	13G24-N240029/900	1
Fagron	Citric acid	14E06-B07293873	1
Fagron	Citric acid	214U-01024	1
Fagron	Citric acid	22031401	1
Fagron	Citric acid	100614d	1
Fagron	Citric acid	13K18-B34-291520	5
Fagron	Citric acid	1714Q1024	1
Fagron	Citric acid	13L18-B39-289347	4
Fagron	Citric acid	13K18-B34-288371	4
Fagron	Citric acid	13K18-B34-291519	5
Fagron	Citric acid	13K18-B34-288372	9
Fagron	Citric acid	13K18-B34-291071	4
Fagron	Citric acid	14g04-b09293877	1
Fagron	Citric acid	700309-0006	1
Fagron	Citric acid	14K05-B04-302816	6
Fagron	Citric acid	13F12-N08	1
Fagron	Citric acid	13K18-B34-293871	2
Fagron	Citric acid	11H30-N05	1
Fagron	Citric acid	13K18-B34-291522	4
Fagron	Clindamycin hydrochlo...	13J03-B02_288654	1
Fagron	Clindamycin hydrochlo...	18D11B01	1
Fagron	Clindamycin hydrochlo...	15E06-B12-310472	2
Fagron	Clindamycin hydrochlo...	14K04-B05-309851	1
Fagron	Clindamycin hydrochlo...	14K04-B05-304996	1
Fagron	Clindamycin hydrochlo...	13G24-N09	8
Fagron	Clindamycin hydrochlo...	13G24N09	1
Fagron	Clindamycin hydrochlo...	13J03-B02-288655	1
Fagron	Clindamycin hydrochlo...	14D23-B03	1
Fagron	Clindamycin hydrochlo...	13J03B02288654	1
Fagron	Clindamycin hydrochlo...	1010616	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clindamycin hydrochlo...	14K04-B05-304995	1
Fagron	Clindamycin hydrochlo...	15L22-B03-318787	1
Fagron	Clindamycin hydrochlo...	2714Q-03534	1
Fagron	Clindamycin hydrochlo...	111105803	1
Fagron	Clindamycin hydrochlo...	15E06-B12-313089	2
Fagron	Clindamycin hydrochlo...	13j03-B02-288653	1
Fagron	Clindamycin hydrochlo...	13C15-N06	1
Fagron	Clindamycin hydrochlo...	13J03-B02-288654	3
Fagron	Clindamycin hydrochlo...	13J03-B02	1
Fagron	Clindamycin hydrochlo...	P003WS13091401	1
Fagron	Clindamycin hydrochlo...	18D11-B01-353162	1
Fagron	Clindamycin hydrochlo...	16D05-B08-327127	6
Fagron	Clindamycin hydrochlo...	15E06-B12-311886	2
Fagron	Clindamycin hydrochlo...	16H03-B02-334344	1
Fagron	Clindamycin hydrochlo...	P003-WS10100301	1
Fagron	Clindamycin hydrochlo...	18041801	1
Fagron	Clindamycin hydrochlo...	16H03-B02-338802	1
Fagron	Clindamycin hydrochlo...	16H03-B02-338803	2
Fagron	Clindamycin hydrochlo...	16H03-B02-334345	1
Fagron	Clindamycin hydrochlo...	16d05-b08-325089	1
Fagron	Clindamycin hydrochlo...	16d05-b08-332107	1
Fagron	Clindamycin hydrochlo...	16D05-B08-332107	1
Fagron	Clindamycin hydrochlo...	16D05-B08-328260	1
Fagron	Clindamycin hydrochlo...	16D05-B08-320327	3
Fagron	Clindamycin hydrochlo...	16D05-B08-325087	1
Fagron	Clindamycin hydrochlo...	16D05-B08-330624	1
Fagron	Clindamycin hydrochlo...	16D05-B08-332105	1
Fagron	Clindamycin hydrochlo...	16D05-B08-323321	1
Fagron	Clindamycin hydrochlo...	16D05-B08	2
Fagron	Clindamycin hydrochlo...	1600324C	1
Fagron	Clindamycin hydrochlo...	18A18-B08	2
Fagron	Clindamycin hydrochlo...	18A18-B08-346107	1
Fagron	Clindamycin hydrochlo...	16H03-B02-336706	2
Fagron	Clindamycin hydrochlo...	14K04-B05-313090	1
Fagron	Clindamycin hydrochlo...	15L22-B03-316670	1
Fagron	Clindamycin hydrochlo...	18D11-B01	2
Fagron	Clindamycin hydrochlo...	18D11-B01-351175	1
Fagron	Clindamycin hydrochlo...	13A22-N01	7
Fagron	Clindamycin hydrochlo...	18A18-B08-347483	3
Fagron	Clindamycin hydrochlo...	15E06-B12-314066	1
Fagron	Clindamycin hydrochlo...	12C19-N06	1
Fagron	Clindamycin hydrochlo...	18D11-B01-351176	1
Fagron	Clindamycin hydrochlo...	16D05-B08-321572	1
Fagron	Clindamycin hydrochlo...	12F19-N07	3
Fagron	Clindamycin hydrochlo...	18D11-B01-349677	1
Fagron	Clindamycin hydrochlo...	18D11-B01-352733	1
Fagron	Clioquinol	14C21-B06-293132	2
Fagron	Clioquinol	14C21-B06-293104	7
Fagron	Clioquinol	13I26-N15	4
Fagron	Clioquinol	2514U2562	1
Fagron	Clioquinol	13DS15-N01	1
Fagron	Clioquinol	13A29-N07-3398310111	1
Fagron	Clioquinol	7D20-N03	1
Fagron	Clioquinol	2014Q2562	1
Fagron	Clioquinol	14052014C	1
Fagron	Clioquinol	13I26N15	1
Fagron	Clioquinol	1205110014	1
Fagron	Clioquinol	28031401	1
Fagron	Clioquinol	14-C21-B06-293132	1
Fagron	Clioquinol	14c21B06293104	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clioquinol	13F12-N20	2
Fagron	Clioquinol	14I11-B03	2
Fagron	Clioquinol	13F12-N200463480	1
Fagron	Clioquinol	1305120028	1
Fagron	Clioquinol	13D15-N01	8
Fagron	Clioquinol	3916U-02562	1
Fagron	Clioquinol	14I11-B03-309181	2
Fagron	Clioquinol	12L04-N08	10
Fagron	Clioquinol	16J18-B04	2
Fagron	Clioquinol	15L30-B05-314974	8
Fagron	Clioquinol	1505120044	1
Fagron	Clioquinol	15L30-B05	2
Fagron	Clioquinol	16/18-B04-328511	1
Fagron	Clioquinol	16J18-B04-329141	3
Fagron	Clioquinol	14/11-Bo3	1
Fagron	Clioquinol	16J18-B04-32914	1
Fagron	Clioquinol	15L30-B05-321899	1
Fagron	Clioquinol	16/18-b04-328511	1
Fagron	Clioquinol	16J18-B04-328520	3
Fagron	Clioquinol	16J18-B04-328511	5
Fagron	Clioquinol	15L30-B05-321899	6
Fagron	Clioquinol	1326-N15	3
Fagron	Clioquinol	15L30-B05-319616	2
Fagron	Clioquinol	14I11-B03-300049	1
Fagron	Clioquinol	15L30-321899	1
Fagron	Clioquinol	12F20-N05	4
Fagron	Clioquinol	1605100008	1
Fagron	Clioquinol	16J18B04	1
Fagron	Clioquinol	1005090024	1
Fagron	Clioquinol	11G20-N09	1
Fagron	Clioquinol	11K21-N18	4
Fagron	Clioquinol	14I11-B03-300049	4
Fagron	Clioquinol	14C21806293104	1
Fagron	Clioquinol	14C21-B06-297159	3
Fagron	Clioquinol	14/11-B03-300049	3
Fagron	Clioquinol	1411-B03-300049	2
Fagron	Clioquinol	1618B04-329141	1
Fagron	Clioquinol	1411B03300049	1
Fagron	Clioquinol	14I11-B03-3000049	1
Fagron	Clioquinol	14I11-B03-309048	2
Fagron	Clioquinol	14/11-B03-309181	4
Fagron	Clioquinol	14/11-803-309181	1
Fagron	Clioquinol	18E08-B06	1
Fagron	Clioquinol	18B05-B07-345483	1
Fagron	Clioquinol	16/18-B04-332852	1
Fagron	Clioquinol	14 11-B03-300049	1
Fagron	Clioquinol	13/26-N15	2
Fagron	Clobetasol propionate	14C03-B07	3
Fagron	Clobetasol propionate	13K08-B30-287727	1
Fagron	Clobetasol propionate	13K08-B30	2
Fagron	Clobetasol propionate	1211021	1
Fagron	Clobetasol propionate	CL0103/0518	1
Fagron	Clobetasol propionate	11129-N07	1
Fagron	Clobetasol propionate	18/12-B02-355065	1
Fagron	Clobetasol propionate	18/12-B02-355390	2
Fagron	Clobetasol propionate	1314705	1
Fagron	Clobetasol propionate	13G23-N01NT0421M	1
Fagron	Clobetasol propionate	13Ao2-N04	1
Fagron	Clobetasol propionate	13G23-N01	2
Fagron	Clobetasol propionate	18I12-B02-353202	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clobetasol propionate	12E08N04	1
Fagron	Clobetasol propionate	3713M-02835	1
Fagron	Clobetasol propionate	NT0421M	3
Fagron	Clobetasol propionate	13G23-N01	17
Fagron	Clobetasol propionate	13A02-N04	18
Fagron	Clobetasol propionate	11H10-N02	1
Fagron	Clobetasol propionate	NS0292-M	1
Fagron	Clobetasol propionate	NS0486M	1
Fagron	Clobetasol propionate	11D26-N04	1
Fagron	Clobetasol propionate	11/29-N07	2
Fagron	Clobetasol propionate	12J25-N01	5
Fagron	Clobetasol propionate	11L12-N05	1
Fagron	Clobetasol propionate	11-05-N06	1
Fagron	Clobetasol propionate	12E08-N04	13
Fagron	Clobetasol propionate	11300703	5
Fagron	Clobetasol propionate	16K24-B04-333143	3
Fagron	Clobetasol propionate	13G23N01	1
Fagron	Clobetasol propionate	16K24-B04-333120	1
Fagron	Clobetasol propionate	18112-B02-353202	1
Fagron	Clobetasol propionate	13K08-B30-287730	2
Fagron	Clobetasol propionate	CL0101/0417	1
Fagron	Clobetasol propionate	17096708	12
Fagron	Clobetasol propionate	17G14-B12-340176	2
Fagron	Clobetasol propionate	17E04-B02-335335	1
Fagron	Clobetasol propionate	16K24B04333121	1
Fagron	Clobetasol propionate	17G14-B12-340175	5
Fagron	Clobetasol propionate	14C03-B07-293562	9
Fagron	Clobetasol propionate	17E04-B02-336882	2
Fagron	Clobetasol propionate	17G14-B12-341742	2
Fagron	Clobetasol propionate	16K24-B04	1
Fagron	Clobetasol propionate	17G14-B12-341743	1
Fagron	Clobetasol propionate	17L18-B02-347992	1
Fagron	Clobetasol propionate	18112-Bo2-353201	1
Fagron	Clobetasol propionate	17L18-B02-347993	4
Fagron	Clobetasol propionate	17E04B02336882	1
Fagron	Clobetasol propionate	13K08-B30-287729	2
Fagron	Clobetasol propionate	17L18-B02-345404	2
Fagron	Clobetasol propionate	17L18B02-344841	1
Fagron	Clobetasol propionate	17L18-B02347991	1
Fagron	Clobetasol propionate	17L18B02347991	1
Fagron	Clobetasol propionate	18I12-B02	1
Fagron	Clobetasol propionate	17L18-B02-345403	1
Fagron	Clobetasol propionate	18112-B02-354809	1
Fagron	Clobetasol propionate	18/12-B02-353202	1
Fagron	Clobetasol propionate	13K08B30287730	1
Fagron	Clobetasol propionate	13k08b30287728	1
Fagron	Clobetasol propionate	17L18-B02-344842	3
Fagron	Clobetasol propionate	16K24-B04-331220	1
Fagron	Clobetasol propionate	17E04-B02-338441	3
Fagron	Clobetasol propionate	16G18-B02-325572	1
Fagron	Clobetasol propionate	15J15-B04-317698	6
Fagron	Clobetasol propionate	15J15-B04-316615	1
Fagron	Clobetasol propionate	15c23-b02-310473	2
Fagron	Clobetasol propionate	1405018-02	1
Fagron	Clobetasol propionate	15C23-B02-310522	2
Fagron	Clobetasol propionate	15H26-B03-312860	1
Fagron	Clobetasol propionate	14G30-B07-302112	2
Fagron	Clobetasol propionate	14C03-B07-305987	2
Fagron	Clobetasol propionate	15a14-b02-303125	2
Fagron	Clobetasol propionate	15C23-B02-310473	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clobetasol propionate	14G30-B07-302112	1
Fagron	Clobetasol propionate	14C03-B07-305985	1
Fagron	Clobetasol propionate	14G30-B07-302113	2
Fagron	Clobetasol propionate	15A14-B02-303125	10
Fagron	Clobetasol propionate	14G30-B07-302114	12
Fagron	Clobetasol propionate	14G30-B07-302116	4
Fagron	Clobetasol propionate	11C16-N01	1
Fagron	Clobetasol propionate	14G30-B07-303114	1
Fagron	Clobetasol propionate	13K08-B03-287728	1
Fagron	Clobetasol propionate	14G30-B07-302115	2
Fagron	Clobetasol propionate	14C03-B07-293561	2
Fagron	Clobetasol propionate	15J15-B04-316616	1
Fagron	Clobetasol propionate	15C23-B02-310473	2
Fagron	Clobetasol propionate	15H26-B03-312355	2
Fagron	Clobetasol propionate	16G18-B02-325571	2
Fagron	Clobetasol propionate	16E04-B12-321034	2
Fagron	Clobetasol propionate	16G18B02325571	1
Fagron	Clobetasol propionate	16E04-B12-321032	1
Fagron	Clobetasol propionate	16G18-B02-32571	1
Fagron	Clobetasol propionate	16H22-B07-330169	3
Fagron	Clobetasol propionate	16G18-B02-325573	7
Fagron	Clobetasol propionate	15j15-B04-317698	1
Fagron	Clobetasol propionate	16E04-B12-321033	2
Fagron	Clobetasol propionate	15/15-B04-317699	1
Fagron	Clobetasol propionate	16G18-B02-324627	6
Fagron	Clobetasol propionate	16H22-B07	1
Fagron	Clobetasol propionate	16E04B12321033	1
Fagron	Clobetasol propionate	16E04-B12-322020	4
Fagron	Clobetasol propionate	16K24-B04-332086	6
Fagron	Clobetasol propionate	16E04-B12-321036	2
Fagron	Clobetasol propionate	15J15B04316616	1
Fagron	Clobetasol propionate	16K24-B04-333120	3
Fagron	Clobetasol propionate	15J15-B04	1
Fagron	Clobetasol-17-propion...	12 19-N01	1
Fagron	Clobetasol-17-propion...	13W16-N01	1
Fagron	Clobetasol-17-propion...	12K08-N02	2
Fagron	Clobetasol-17-propion...	12I19-N01	1
Fagron	Clobetasol-17-propion...	12/19-N01003637	1
Fagron	Clobetasol-17-propion...	12L19-N01	1
Fagron	Clobetasol-17-propion...	12/19-No1	1
Fagron	Clobetasol-17-propion...	12L19-N06	1
Fagron	Clonidine hydrochloride	14H11-B04-298642	1
Fagron	Clonidine hydrochloride	14H11-B04-311163	42
Fagron	Clotrimazole	12E30-N02	1
Fagron	Clotrimazole	14b05-b35-294935	1
Fagron	Clotrimazole	15G14-B06-311055	1
Fagron	Clotrimazole	14b05-b35-299941	1
Fagron	Clotrimazole	14C04-B312-299959	1
Fagron	Clotrimazole	14C04-B31-299959	2
Fagron	Clotrimazole	14B06-B06	1
Fagron	Clotrimazole	15G14-B05-311100	10
Fagron	Clotrimazole	16B09-B06-311785	4
Fagron	Clotrimazole	15G14B06311054	1
Fagron	Clotrimazole	15G14-B06-311054	4
Fagron	Clotrimazole	14B05-B35-299944	2
Fagron	Clotrimazole	14C04B31	1
Fagron	Clotrimazole	11A05-N14	2
Fagron	Clotrimazole	14B05-B36-291166	3
Fagron	Clotrimazole	14B05-B35-299941	2
Fagron	Clotrimazole	14B05-B36-291167	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clotrimazole	CLT/1307/0065A	2
Fagron	Clotrimazole	14b05-B36-291167	1
Fagron	Clotrimazole	14B06-B06-299942	1
Fagron	Clotrimazole	14C04-B32-291657	1
Fagron	Clotrimazole	14b05b36291166	1
Fagron	Clotrimazole	14223201	1
Fagron	Clotrimazole	14B05-B35-293543	1
Fagron	Clotrimazole	14C04-B31-294939	7
Fagron	Clotrimazole	14B05-B35-293544	5
Fagron	Clotrimazole	14C04-B31-291660	2
Fagron	Clotrimazole	14C04-B30-291656	11
Fagron	Clotrimazole	14C04-B31	6
Fagron	Clotrimazole	14b05-b36-291166	1
Fagron	Clotrimazole	14B05-B35	3
Fagron	Clotrimazole	14B05-B36-291165	2
Fagron	Clotrimazole	14B05-B36-2991166	1
Fagron	Clotrimazole	14B05B36291167	1
Fagron	Clotrimazole	14C04-B31-294951	1
Fagron	Clotrimazole	14B05-B35-294932	1
Fagron	Clotrimazole	14C04-B31-299954	4
Fagron	Clotrimazole	14b05-b35-293544	1
Fagron	Clotrimazole	14B05-B35-294935	4
Fagron	Clotrimazole	14C04-B31-294940	1
Fagron	Clotrimazole	14C04-B31-299960	6
Fagron	Clotrimazole	14E26-B02	1
Fagron	Clotrimazole	14B05-B35-293542	3
Fagron	Clotrimazole	14B05B35294935	1
Fagron	Clotrimazole	14B05-B36-291168	3
Fagron	Clotrimazole	14B05-B36	8
Fagron	Clotrimazole	11K21-N19	2
Fagron	Clotrimazole	12A05-N10	3
Fagron	Clotrimazole	12E02-N12	10
Fagron	Clotrimazole	11L12-N10	4
Fagron	Clotrimazole	12a05-n10	1
Fagron	Clotrimazole	20111018	4
Fagron	Clotrimazole	12I21N04	1
Fagron	Clotrimazole	12I21-N04	8
Fagron	Clotrimazole	20151112-3	1
Fagron	Clotrimazole	12e02-N12	1
Fagron	Clotrimazole	12/21-N04	3
Fagron	Clotrimazole	C22H17CIN2	1
Fagron	Clotrimazole	11G15-N10	1
Fagron	Clotrimazole	14C04B32291657	1
Fagron	Clotrimazole	12I21-N04	1
Fagron	Clotrimazole	12E30-N07	13
Fagron	Clotrimazole	20120603	3
Fagron	Clotrimazole	12/21-N01	6
Fagron	Clotrimazole	12 21-N01	2
Fagron	Clotrimazole	1221n04	1
Fagron	Clotrimazole	14B06-B06-299943	6
Fagron	Clotrimazole	12I21-N01	7
Fagron	Clotrimazole	9041304	1
Fagron	Clotrimazole	4812E-01276 (NP)	1
Fagron	Clotrimazole	13C25-N10	10
Fagron	Clotrimazole	1221N04	1
Fagron	Clotrimazole	12I16-N04	1
Fagron	Clotrimazole	12K15-N01	1
Fagron	Clotrimazole	12F02-N12	1
Fagron	Clotrimazole	15G14-B05-311103	8
Fagron	Clotrimazole	11E17-N10	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clotrimazole	15G14-Bo5-311101	1
Fagron	Clotrimazole	17A19-B04-336568	1
Fagron	Clotrimazole	16J07-B02-333606	2
Fagron	Clotrimazole	17J25-B06-345407	1
Fagron	Clotrimazole	17J25-B06-345045	1
Fagron	Clotrimazole	16J07-B02-341842	1
Fagron	Clotrimazole	16607-B02-333606	1
Fagron	Clotrimazole	16J07-B02-344816	1
Fagron	Clotrimazole	17J25-B06-345044	2
Fagron	Clotrimazole	17360201	17
Fagron	Clotrimazole	16J07-B02-344818	5
Fagron	Clotrimazole	16J07-B02-340513	3
Fagron	Clotrimazole	16J07-B02-340513-201	1
Fagron	Clotrimazole	17A19-B04-340476	2
Fagron	Clotrimazole	17A19-B04-338577	2
Fagron	Clotrimazole	16J07-B02-340511	2
Fagron	Clotrimazole	17A19-B04-334768	1
Fagron	Clotrimazole	17A19-SB04-332097	1
Fagron	Clotrimazole	16J07-B02-338283	4
Fagron	Clotrimazole	16J07-B02-337573	3
Fagron	Clotrimazole	CLT/150570016b	1
Fagron	Clotrimazole	17J25-B06	2
Fagron	Clotrimazole	18A18-B05-344819	1
Fagron	Clotrimazole	16J07-B02-347864	4
Fagron	Clotrimazole	18A18-B05-347865	4
Fagron	Clotrimazole	18G10-B06-352708	1
Fagron	Clotrimazole	18A18-B05-350781	1
Fagron	Clotrimazole	18A18-B05-352114	2
Fagron	Clotrimazole	18A18-B05-352115	2
Fagron	Clotrimazole	18A18B05	1
Fagron	Clotrimazole	16B09-B06-311507	5
Fagron	Clotrimazole	18A18-B05-350780	2
Fagron	Clotrimazole	19B14-B01-359532	1
Fagron	Clotrimazole	16J07B02341842	2
Fagron	Clotrimazole	18A18B05352114	1
Fagron	Clotrimazole	700135-0003	1
Fagron	Clotrimazole	18G10-B06-353348	1
Fagron	Clotrimazole	18A18-B05-350778	4
Fagron	Clotrimazole	700120-0003	1
Fagron	Clotrimazole	18a18-b05-348429	1
Fagron	Clotrimazole	1A1-B05-34765	1
Fagron	Clotrimazole	17J25-B06-349161	2
Fagron	Clotrimazole	CLT/1711/0077A	1
Fagron	Clotrimazole	18A18-B05-348429	2
Fagron	Clotrimazole	17J25-B06-347476	2
Fagron	Clotrimazole	16J07-B02-344817	3
Fagron	Clotrimazole	20171127Ma1	1
Fagron	Clotrimazole	16234521	8
Fagron	Clotrimazole	18G10-B06352708	1
Fagron	Clotrimazole	16J07-B02-329386	4
Fagron	Clotrimazole	20170427Ma9	1
Fagron	Clotrimazole	16B11-B04-327611	3
Fagron	Clotrimazole	CLT/1604/00227A	1
Fagron	Clotrimazole	CLT/160/00058	1
Fagron	Clotrimazole	16E20-B10-326010	3
Fagron	Clotrimazole	16E20-B10-327349	6
Fagron	Clotrimazole	160704602	3
Fagron	Clotrimazole	16E20-B10-321597	2
Fagron	Clotrimazole	15G14-B05-311102	3
Fagron	Clotrimazole	16E20-B10-321598	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Clotrimazole	16E20-B10-32197	1
Fagron	Clotrimazole	15G14-B06-31105	1
Fagron	Clotrimazole	16B09-B06	1
Fagron	Clotrimazole	17A19-B04332114	1
Fagron	Clotrimazole	3716Q-01276	1
Fagron	Clotrimazole	16J07-B02-330216	7
Fagron	Clotrimazole	16B09-B06-311506	4
Fagron	Clotrimazole	2616A-01276	1
Fagron	Clotrimazole	15G14-B06-311056	1
Fagron	Clotrimazole	15G14-B05	5
Fagron	Clotrimazole	17A19-B04_332114	1
Fagron	Clotrimazole	CLT/1604/0027A	2
Fagron	Colistin sulfate	15E15-B01-315340	1
Fagron	Colistin sulfate	16F03-B01	1
Fagron	Colophony	15186807	8
Fagron	Copper (II) sulfate p...	13A25-N06	1
Fagron	Corn starch	15310902	2
Fagron	Cyanocobalamin	13124-B01-292062	1
Fagron	Cyanocobalamin	15H31-B02-315113	1
Fagron	Cyanocobalamin	17A27-B04-339672	1
Fagron	Cyanocobalamin	17A27-B04	1
Fagron	Cyanocobalamin	10K02-N05	1
Fagron	Cyanocobalamin	10K02-N05	2
Fagron	Cyproterone acetate	8B11-N06	1
Fagron	Cyproterone acetate	11L22-N02	1
Fagron	Cyproterone acetate	12F25-N09	3
Fagron	Cyproterone acetate	12D02-N13	1
Fagron	Cyproterone acetate	B0392/110020	3
Fagron	Cyproterone acetate	12J24-N02	1
Fagron	Cyproterone acetate	18G25-B02-352307	1
Fagron	Cyproterone acetate	18G25-B02-352728	3
Fagron	Cyproterone acetate	13/16-N04	2
Fagron	Cyproterone acetate	17D11-B07	1
Fagron	Cyproterone acetate	17D11-B07-336729	2
Fagron	Cyproterone acetate	17D11-B07-338440	1
Fagron	Cyproterone acetate	EU-CPA-170301	1
Fagron	Cyproterone acetate	16K09-B04-329037	1
Fagron	Cyproterone acetate	15I11-B01-315381	2
Fagron	Cyproterone acetate	13j16-no4	1
Fagron	Cyproterone acetate	13j16n04	1
Fagron	Cyproterone acetate	13J16-N04	2
Fagron	Cyproterone acetate	13/09-N03	1
Fagron	Cyproterone acetate	12L19-N05	3
Fagron	Cyproterone acetate	18G25-B02-355388	1
Fagron	Cystine	15F08-B05-310685	1
Fagron	Cystine	12D26-B04-285823	2
Fagron	Deoxycholic acid	13J21-N14	1
Fagron	Dequalinium chloride	1407029-01	1
Fagron	Dequalinium chloride	12G05-N07	3
Fagron	Dequalinium chloride	12D05-N09	1
Fagron	Dequalinium chloride	15B16-B04-346969	3
Fagron	Dequalinium chloride	13A23-N02	3
Fagron	Dequalinium chloride	15B16-B04-340177	2
Fagron	Dequalinium chloride	9245-90.81	1
Fagron	Dequalinium chloride	12B14-N01	1
Fagron	Dequalinium chloride	15B16-B04-316267	2
Fagron	Dequalinium chloride	15B16-B04	1
Fagron	Dequalinium chloride	15B16-B04-304432	1
Fagron	Dequalinium chloride	13I26-N10	1
Fagron	Dequalinium chloride	14C10-B03-292286	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Dequalinium chloride	13D26-N06	3
Fagron	Dequalinium chloride	14C10-B03	1
Fagron	Dequalinium chloride	15B16B04304432	1
Fagron	Dequalinium chloride	10H05-N07	1
Fagron	Dexamethasone	15I04-B04-315528	1
Fagron	Dexamethasone	16I06-B03-327488	3
Fagron	Dexamethasone	14A22-B03-290701	1
Fagron	Dexamethasone	13D12-N09	4
Fagron	Dexamethasone	2314A-01327	1
Fagron	Dexamethasone	15A05-B01-034636	1
Fagron	Dexamethasone	700107-0002	1
Fagron	Dexamethasone	1209008-01	11
Fagron	Dexamethasone	12B06-N12	2
Fagron	Dexamethasone	12H21-N11	2
Fagron	Dexamethasone	11L08-N01	4
Fagron	Dexamethasone	12D26-N09	3
Fagron	Dexamethasone	16K24B05330281	1
Fagron	Dexamethasone	15A05-B01-304636	3
Fagron	Dexamethasone	16/06-B03-327323	1
Fagron	Dexamethasone	16I06-B03-327323	1
Fagron	Dexamethasone	15K17-B07-321457	1
Fagron	Dexamethasone	13L30-B34-289691	1
Fagron	Dexamethasone	15C16-B07	1
Fagron	Dexamethasone	15/04-B04-315528	1
Fagron	Dexamethasone	14D11-B01-294032	1
Fagron	Dexamethasone	15A05-B01	1
Fagron	Dexamethasone	11F30-N02	1
Fagron	Dexamethasone	14A22-B03	2
Fagron	Dexamethasone	14A22-B03-294030	1
Fagron	Dexamethasone	17C13-B07-332461	1
Fagron	Dexamethasone	17L19-B03-344396	1
Fagron	Dexamethasone	11D26-N10	3
Fagron	Dexamethasone	18D09-B06-348185	1
Fagron	Dexamethasone	15C16-B07-308252	1
Fagron	Dexamethasone	17L19-B03-344100	2
Fagron	Dexamethasone	NUD160702	1
Fagron	Dexamethasone	17E31-B01-331860	1
Fagron	Dexamethasone	DZ11005K12B06-N12	1
Fagron	Dexamethasone	13G08-N13	5
Fagron	Dexamethasone	17104203	5
Fagron	Dexamethasone	13D12-N19	1
Fagron	Dexamethasone	17E31-B01-338016	1
Fagron	Dexamethasone	16K24-B05-330281	2
Fagron	Dexamethasone 1:10 (r...	12L19-N06	1
Fagron	Dexamethasone 1:10 (r...	12K06-N09	1
Fagron	Dexamethasone 1:10 (r...	12C20-N03	4
Fagron	Dexamethasone acetate	1110007	1
Fagron	Dextromethorphan hydr...	16A22-B03-321371	5
Fagron	Dextromethorphan hydr...	14H27-B04-299502	7
Fagron	Dextromethorphan hydr...	16A22-B03-319331	1
Fagron	D-galactose	16G22-B01-329906	1
Fagron	Diclofenac sodium	17J19-B01-350960	3
Fagron	Diclofenac sodium	17F02-B05	1
Fagron	Diclofenac sodium	17F02	1
Fagron	Diclofenac sodium	16129	1
Fagron	Diclofenac sodium	17J19-B01	1
Fagron	Diclofenac sodium	18H24-B07-357416	1
Fagron	Diclofenac sodium	17j19-B01-342310	1
Fagron	Diclofenac sodium	13E03-N02	1
Fagron	Diclofenac sodium	5041401	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Diclofenac sodium	15K23-B03-327802	5
Fagron	Diclofenac sodium	15K23-B03-325753	4
Fagron	Diclofenac sodium	16I29-B04-332372	2
Fagron	Diclofenac sodium	15K23-B03-321673	3
Fagron	Diclofenac sodium	15K23-B03-321673	1
Fagron	Diclofenac sodium	17F02-B05-338797	6
Fagron	Diclofenac sodium	15K23-B03	2
Fagron	Diclofenac sodium	15K23-B03-318519	4
Fagron	Diclofenac sodium	15D14-B04-315835	1
Fagron	Diclofenac sodium	15D14-B04-314450	3
Fagron	Diclofenac sodium	14031411	1
Fagron	Diclofenac sodium	15K23-B03-324025	3
Fagron	Diclofenac sodium	14G17B08297032	1
Fagron	Diclofenac sodium	12L28-B05-289773	2
Fagron	Diclofenac sodium	RDS/1405/0159A	1
Fagron	Diclofenac sodium	15D14-B04	1
Fagron	Diltiazem hydrochloride	12I12-N02	4
Fagron	Diltiazem hydrochloride	11E18-N04	3
Fagron	Diltiazem hydrochloride	12K07-N02	2
Fagron	Diltiazem hydrochloride	12D05-N070	1
Fagron	Diltiazem hydrochloride	12d05-n07	1
Fagron	Diltiazem hydrochloride	12/12-N02	1
Fagron	Diltiazem hydrochloride	11I01-N08	1
Fagron	Diltiazem hydrochloride	10D28-N23	1
Fagron	Diltiazem hydrochloride	11-01-N08	1
Fagron	Diltiazem hydrochloride	612213562	1
Fagron	Diltiazem hydrochloride	12F19-N10	4
Fagron	Diltiazem hydrochloride	12D05-N07	13
Fagron	Diltiazem hydrochloride	12/12-N02	8
Fagron	Diltiazem hydrochloride	3713M-03562	1
Fagron	Diltiazem hydrochloride	10808-N01	1
Fagron	Diltiazem hydrochloride	17D11-B11-333704	12
Fagron	Diltiazem hydrochloride	16G20-B03-328253	2
Fagron	Diltiazem hydrochloride	16L19-B01-330396	7
Fagron	Diltiazem hydrochloride	AZ-0100716	1
Fagron	Diltiazem hydrochloride	16L19-B01	1
Fagron	Diltiazem hydrochloride	17D-11-B11-333704	1
Fagron	Diltiazem hydrochloride	DIL/M-21816	4
Fagron	Diltiazem hydrochloride	16GH20-B03-325903	1
Fagron	Diltiazem hydrochloride	16c14-b11-322152	1
Fagron	Diltiazem hydrochloride	17D11-B11-334837	13
Fagron	Diltiazem hydrochloride	17D11-B11-336727	3
Fagron	Diltiazem hydrochloride	17D11-B11-334192	2
Fagron	Diltiazem hydrochloride	17D11-B11-333705	5
Fagron	Diltiazem hydrochloride	16L19-B01_330396	1
Fagron	Diltiazem hydrochloride	17D11-B11	5
Fagron	Diltiazem hydrochloride	16L19-B01-330365	11
Fagron	Diltiazem hydrochloride	16g20-b03-325903	1
Fagron	Diltiazem hydrochloride	16C14B11322152	1
Fagron	Diltiazem hydrochloride	16C14-B11-320568	5
Fagron	Diltiazem hydrochloride	16C14B11320849	1
Fagron	Diltiazem hydrochloride	16a18-b07-318732	2
Fagron	Diltiazem hydrochloride	16C14-B11-322642	16
Fagron	Diltiazem hydrochloride	16C14-B11-322153	8
Fagron	Diltiazem hydrochloride	16C14B11320568	2
Fagron	Diltiazem hydrochloride	16c14-b11-211642	1
Fagron	Diltiazem hydrochloride	16C14-B11-322152	5
Fagron	Diltiazem hydrochloride	16A18B07316668	1
Fagron	Diltiazem hydrochloride	16C14-B11	3
Fagron	Diltiazem hydrochloride	16G20B03325903	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Diltiazem hydrochloride	16G20-B03	4
Fagron	Diltiazem hydrochloride	16G20-B03-328133	13
Fagron	Diltiazem hydrochloride	13A02-N05	10
Fagron	Diltiazem hydrochloride	16A14-B11-322153	1
Fagron	Diltiazem hydrochloride	17H25-B01-340697	13
Fagron	Diltiazem hydrochloride	16L19-B01-330395	1
Fagron	Diltiazem hydrochloride	17H25-B01-340699	2
Fagron	Diltiazem hydrochloride	AZ-0641216	1
Fagron	Diltiazem hydrochloride	16G20B03-325903	1
Fagron	Diltiazem hydrochloride	18A29B01347886	1
Fagron	Diltiazem hydrochloride	18A29-B01-350964	1
Fagron	Diltiazem hydrochloride	18A29-B01-345623	1
Fagron	Diltiazem hydrochloride	16L19B01330396	1
Fagron	Diltiazem hydrochloride	18A29-B01-352909	1
Fagron	Diltiazem hydrochloride	18A29-B01-358861	1
Fagron	Diltiazem hydrochloride	18A29-B01-352723	1
Fagron	Diltiazem hydrochloride	18A29-B01-352907	2
Fagron	Diltiazem hydrochloride	18A29-B01-32907	1
Fagron	Diltiazem hydrochloride	18H09-B03-354798	3
Fagron	Diltiazem hydrochloride	18H09B03357184	1
Fagron	Diltiazem hydrochloride	18A-29-B01	1
Fagron	Diltiazem hydrochloride	18H09-B03	2
Fagron	Diltiazem hydrochloride	18H09-B03-357184	1
Fagron	Diltiazem hydrochloride	18A29-B01-348861	14
Fagron	Diltiazem hydrochloride	16C14-B11-320849	15
Fagron	Diltiazem hydrochloride	17H25-B01-341598	1
Fagron	Diltiazem hydrochloride	18A29-B01-347886	8
Fagron	Diltiazem hydrochloride	17H25-B01-340697	1
Fagron	Diltiazem hydrochloride	17d11-b11-336727	1
Fagron	Diltiazem hydrochloride	17H25-B01-341597	2
Fagron	Diltiazem hydrochloride	17H25-B01	2
Fagron	Diltiazem hydrochloride	17H25-B01-341596	8
Fagron	Diltiazem hydrochloride	17H25B01-341598	1
Fagron	Diltiazem hydrochloride	17H25-B01-343273	5
Fagron	Diltiazem hydrochloride	18A29-B01-345623	4
Fagron	Diltiazem hydrochloride	18A29-B01-345625	1
Fagron	Diltiazem hydrochloride	17H25-B01-340698	1
Fagron	Diltiazem hydrochloride	18A29-B01	5
Fagron	Diltiazem hydrochloride	18a29b01345623	1
Fagron	Diltiazem hydrochloride	18A29-B01-346695	2
Fagron	Diltiazem hydrochloride	18A29-B01-346694	5
Fagron	Diltiazem hydrochloride	17H25-B01-343272	1
Fagron	Diltiazem hydrochloride	18A29-01-347886	1
Fagron	Diltiazem hydrochloride	16A18-B07-318732	18
Fagron	Diltiazem hydrochloride	DIL/M-20215	1
Fagron	Diltiazem hydrochloride	16A18B07316667	1
Fagron	Diltiazem hydrochloride	135968-0009	2
Fagron	Diltiazem hydrochloride	14C12-B04-294134	1
Fagron	Diltiazem hydrochloride	14c12-b04-294134	1
Fagron	Diltiazem hydrochloride	16A18-B07-316667	8
Fagron	Diltiazem hydrochloride	14C12-B04-294133	5
Fagron	Diltiazem hydrochloride	14C12_B04-294133	1
Fagron	Diltiazem hydrochloride	DIL/M-01314	1
Fagron	Diltiazem hydrochloride	14c12-b04-292759	1
Fagron	Diltiazem hydrochloride	14C12B04	1
Fagron	Diltiazem hydrochloride	DIL/M-0314	1
Fagron	Diltiazem hydrochloride	22082014A	1
Fagron	Diltiazem hydrochloride	14C12-B04-292759	5
Fagron	Diltiazem hydrochloride	13346915	9
Fagron	Diltiazem hydrochloride	14C12_B04_292757	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Diltiazem hydrochloride	14C12-B04	7
Fagron	Diltiazem hydrochloride	13J24-B04-286570	1
Fagron	Diltiazem hydrochloride	13J24-B04	3
Fagron	Diltiazem hydrochloride	11J28-N05	1
Fagron	Diltiazem hydrochloride	13J24-B04-286571	2
Fagron	Diltiazem hydrochloride	13J24-B04-286569	12
Fagron	Diltiazem hydrochloride	9J14-N17	1
Fagron	Diltiazem hydrochloride	13J24B04286569	3
Fagron	Diltiazem hydrochloride	12112-N02	1
Fagron	Diltiazem hydrochloride	13G16-N05	24
Fagron	Diltiazem hydrochloride	13A02N05	3
Fagron	Diltiazem hydrochloride	13a02-n05	1
Fagron	Diltiazem hydrochloride	12F19-N100102/574	1
Fagron	Diltiazem hydrochloride	13F13-N05	11
Fagron	Diltiazem hydrochloride	14C12-B04-292757	24
Fagron	Diltiazem hydrochloride	14G22-B08-299781	14
Fagron	Diltiazem hydrochloride	16G20-B03-325903	19
Fagron	Diltiazem hydrochloride	14G22-B08-299783	1
Fagron	Diltiazem hydrochloride	15B12-B01-305399	24
Fagron	Diltiazem hydrochloride	15g20-b08-314375	1
Fagron	Diltiazem hydrochloride	15b12-b01-305399	2
Fagron	Diltiazem hydrochloride	15B12-B01_305399	1
Fagron	Diltiazem hydrochloride	15B12-B01	8
Fagron	Diltiazem hydrochloride	15B12B01305399	2
Fagron	Diltiazem hydrochloride	15G20-B08	3
Fagron	Diltiazem hydrochloride	14G22B08	1
Fagron	Diltiazem hydrochloride	DIL/M-13114	1
Fagron	Diltiazem hydrochloride	15G20-B08-312065	15
Fagron	Diltiazem hydrochloride	15G20-B08-310792	1
Fagron	Diltiazem hydrochloride	15B12-b01	1
Fagron	Diltiazem hydrochloride	15B12-B01-305401	3
Fagron	Diltiazem hydrochloride	15B12-B01-305299	1
Fagron	Diltiazem hydrochloride	15b12-b01-305400	2
Fagron	Diltiazem hydrochloride	15B12/B01305399	1
Fagron	Diltiazem hydrochloride	14G22-B08	3
Fagron	Diltiazem hydrochloride	14G22-B08-299782	6
Fagron	Diltiazem hydrochloride	16A18-B07	5
Fagron	Diltiazem hydrochloride	15G20-B08-314375	8
Fagron	Diltiazem hydrochloride	15B12B01305400	1
Fagron	Diltiazem hydrochloride	15B12-B01-305400	8
Fagron	Diltiazem hydrochloride	DIL/M-10515	2
Fagron	Diltiazem hydrochloride	14G22-B08299781	2
Fagron	Diltiazem hydrochloride	15G20-B02-310792	1
Fagron	Diphenhydramine hydro...	16C25-B03-331026	1
Fagron	Diphenhydramine hydro...	12F28-N01	1
Fagron	Diphenhydramine hydro...	12J04-N02	1
Fagron	Diphenhydramine hydro...	13A14-B40-293346	1
Fagron	Diphenylcyclopropenone	35-16-I-95P3117	1
Fagron	Diphenylcyclopropenone	16F02B10324263	1
Fagron	Diphenylcyclopropenone	1(17B13-B02-34977)	1
Fagron	Diphenylcyclopropenone	18I03B04	1
Fagron	Diphenylcyclopropenone	16F02-B10-327812	3
Fagron	Dithranol	11L20-N13	1
Fagron	Dithranol	17710I-01274	1
Fagron	Dithranol	13A30-N01	4
Fagron	Doxycycline hyclate	17112201	1
Fagron	Doxycycline hyclate	18D18-B01-349567	1
Fagron	Doxycycline hyclate	14K28-B05-312010	1
Fagron	Doxycycline hyclate	16I06-B02	1
Fagron	Doxycycline hyclate	16I06-B02-330167	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Doxycycline hyclate	1603008	1
Fagron	Doxycycline hyclate	20101701	1
Fagron	Doxycycline hyclate	16106-Bo2	1
Fagron	Doxycycline hyclate	16106-B02-330167	1
Fagron	Doxycycline hyclate	16E04-B07-231917	1
Fagron	Doxycycline hyclate	12G12-N15	1
Fagron	Doxycycline hyclate	13D24-B40-298948	1
Fagron	Doxylamine hydrogen s...	11G20-N10	1
Fagron	Doxylamine hydrogen s...	12811420001	1
Fagron	Doxylamine hydrogen s...	12E04-N04	3
Fagron	Doxylamine hydrogen s...	12811328003	1
Fagron	Doxylamine hydrogen s...	13/26-N14	1
Fagron	Doxylamine hydrogen s...	14F03-B01-296845	1
Fagron	Doxylamine hydrogen s...	13E29-N03	4
Fagron	Doxylamine hydrogen s...	12K21-N01	2
Fagron	Doxylamine hydrogen s...	11K17-N04	3
Fagron	Doxylamine hydrogen s...	12G10-N03	6
Fagron	Doxylamine hydrogen s...	14K27-B07-303356	1
Fagron	Doxylamine hydrogen s...	15K04-B01-317386	1
Fagron	Doxylamine hydrogen s...	14F03B01293498	1
Fagron	Doxylamine hydrogen s...	16F06-B02-357986	2
Fagron	Doxylamine hydrogen s...	16F06-B02-353931	4
Fagron	Doxylamine hydrogen s...	18041818	1
Fagron	Doxylamine hydrogen s...	15F01-B04-308354	1
Fagron	Doxylamine hydrogen s...	15K04-B01	1
Fagron	Doxylamine hydrogen s...	15K04-B01-327356	1
Fagron	Doxylamine hydrogen s...	15K04-B01-318794	5
Fagron	Doxylamine hydrogen s...	14K27-B07	1
Fagron	Doxylamine hydrogen s...	483875887	1
Fagron	Doxylamine hydrogen s...	15K04-B01-322165	7
Fagron	Doxylamine hydrogen s...	16F06-B-02-338804	1
Fagron	Drofenine hydrochloride	12D04-N03	1
Fagron	Drofenine hydrochloride	11J24-N01	1
Fagron	Econazole nitrate	13D22-N06	1
Fagron	Econazole nitrate	11F09-N16	1
Fagron	Econazole nitrate	14G22-B07-306926	1
Fagron	Emulsifying cetostear...	12G02-N04	6
Fagron	Emulsifying cetostear...	18A23-B04-352987	1
Fagron	Emulsifying cetostear...	11499002	1
Fagron	Emulsifying cetostear...	14G31-B03-305612	4
Fagron	Emulsifying cetostear...	14G31-B03	1
Fagron	Emulsifying cetostear...	13L17-B05-294999	1
Fagron	Emulsifying cetostear...	13L17-B05-289805	1
Fagron	Emulsifying cetostear...	13L17-B05-294999	2
Fagron	Emulsifying cetostear...	3013M-01078	1
Fagron	Emulsifying cetostear...	13L17B05289805	1
Fagron	Emulsifying cetostear...	13L17-B05-289806	2
Fagron	Emulsifying cetostear...	9652865	1
Fagron	Emulsifying cetostear...	13I06-N04	2
Fagron	Emulsifying cetostear...	13L17-B05-289805	6
Fagron	Emulsifying cetostear...	13F07-N05	12
Fagron	Emulsifying cetostear...	12G02N04	1
Fagron	Emulsifying cetostear...	12602-N04	1
Fagron	Emulsifying cetostear...	8256792	2
Fagron	Emulsifying cetostear...	20260712-2	1
Fagron	Emulsifying cetostear...	11L22-N09	1
Fagron	Emulsifying cetostear...	17F21-B01-344775	2
Fagron	Emulsifying cetostear...	15C31-B99-307062	5
Fagron	Emulsifying cetostear...	15D24-B02-308261	3
Fagron	Emulsifying cetostear...	13L17-B05	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Emulsifying cetostear...	15D24-B02-311993	6
Fagron	Emulsifying cetostear...	18003401	1
Fagron	Emulsifying cetostear...	15C31-B00-307062	1
Fagron	Emulsifying cetostear...	17D28-B06-338141	4
Fagron	Emulsifying cetostear...	17D28-B06	1
Fagron	Emulsifying cetostear...	171216CAE/0016830715	1
Fagron	Emulsifying cetostear...	17D28-B06-334609	3
Fagron	Emulsifying cetostear...	16/13-B01-327355	2
Fagron	Emulsifying cetostear...	16I13-B01	2
Fagron	Emulsifying cetostear...	15789271	2
Fagron	Emulsifying cetostear...	16/13-B01-327864	1
Fagron	Emulsifying cetostear...	17F21-B01-341021	5
Fagron	Emulsifying cetostear...	16B23-B07	1
Fagron	Emulsifying cetostear...	16B23-B07-323025	3
Fagron	Emulsifying cetostear...	16B23-B07-318218	2
Fagron	Emulsifying cetostear...	161006CSA	1
Fagron	Emulsifying cetostear...	16B23-B07-320384	5
Fagron	Emulsifying cetostear...	2215I-01078	1
Fagron	Emulsifying cetostear...	15D24-B02	1
Fagron	Emulsifying cetostear...	16I13-B01_326511	1
Fagron	Emulsifying cetostear...	15D24-B02-308261	1
Fagron	Eosin Y	13fF14-T01-017395	1
Fagron	Eosin Y	17372-87-1	1
Fagron	Eosin Y	17K21-T07-042542	2
Fagron	Eosin Y	15G01-T03-034438	1
Fagron	Eosin Y	15G01-T03-025958	3
Fagron	Eosin Y	12A02-N09	1
Fagron	Eosin Y	15G01-T03-027933	2
Fagron	Eosin Y	11F21-N10	1
Fagron	Eosin Y	13/11-N03	1
Fagron	Eosin Y	EYWS/16/11	2
Fagron	Eosin Y	12H28-N08	2
Fagron	Eosin Y	13/11-N03	3
Fagron	Eosin Y	EYWS/15/13	1
Fagron	Eosin Y	14H28-M09	2
Fagron	Eosin Y	12A02-N09	1
Fagron	Eosin Y	15G01T03025959	1
Fagron	Eosin Y	11E25-N01	1
Fagron	Eosin Y	17K21-T07-062160	1
Fagron	Eosin Y	12H28-N08	1
Fagron	Eosin Y	54710	1
Fagron	Eosin Y	EYWS/09/17	1
Fagron	Eosin Y	17K21-T07-057542	1
Fagron	Eosin Y	17K21T07-054710	1
Fagron	Ephedrine hydrochloride	13L30-B36	1
Fagron	Ephedrine hydrochloride	12017220	1
Fagron	Ephedrine hydrochloride	13B19-N04	1
Fagron	Epinephrine hydrogen ...	15A13-B02-305559	2
Fagron	Epinephrine hydrogen ...	13I04-N03	1
Fagron	Epinephrine hydrogen ...	14E12-B02	1
Fagron	Epinephrine hydrogen ...	12G11-No70076/130	1
Fagron	Epinephrine hydrogen ...	12011902	1
Fagron	Epinephrine hydrogen ...	12G11-N07	1
Fagron	Epinephrine hydrogen ...	16e09-b11-319882	1
Fagron	Epinephrine hydrogen ...	16E09-B11-319882	4
Fagron	Epinephrine hydrogen ...	12/02-N03	1
Fagron	Epinephrine hydrogen ...	13I04-N03	1
Fagron	Erythromycin	14I30-B01-311166	2
Fagron	Erythromycin	1430-b01-300060	1
Fagron	Erythromycin	14I30-B01-308789	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Erythromycin	16G06-B03-329157	4
Fagron	Erythromycin	14I30-B01-308188	4
Fagron	Erythromycin	14/30-B01_308790	1
Fagron	Erythromycin	1430-B01-308789	1
Fagron	Erythromycin	14I30-B01-308189	1
Fagron	Erythromycin	14I30-B01-306927	3
Fagron	Erythromycin	300059	1
Fagron	Erythromycin	EB06814	1
Fagron	Erythromycin	14/30-B01-308789	3
Fagron	Erythromycin	3000-50	1
Fagron	Erythromycin	EB068/14	1
Fagron	Erythromycin	14I30-B01-300050	1
Fagron	Erythromycin	1430B01300059	1
Fagron	Erythromycin	14 30-B01-311166	1
Fagron	Erythromycin	14i30-b01-311166	1
Fagron	Erythromycin	14I30B01308189	1
Fagron	Erythromycin	16/24-B10-329158	1
Fagron	Erythromycin	16<g06-B03-324464	1
Fagron	Erythromycin	16G06-B03324464	1
Fagron	Erythromycin	16G06-B03-329149	4
Fagron	Erythromycin	16G06-B03-324462	2
Fagron	Erythromycin	16G06-B03_0003	1
Fagron	Erythromycin	16g06-b03-324463	2
Fagron	Erythromycin	16G06-B03	7
Fagron	Erythromycin	16G06-B03-324465	6
Fagron	Erythromycin	OEBA/0010516	2
Fagron	Erythromycin	16G06-B03-324463	1
Fagron	Erythromycin	16G06-B03-324464	13
Fagron	Erythromycin	14/30-B01	1
Fagron	Erythromycin	16J24-B10-329154	1
Fagron	Erythromycin	16g06-b03-329150	1
Fagron	Erythromycin	16G06-B03-324463	10
Fagron	Erythromycin	16G06B03-329149	1
Fagron	Erythromycin	16J24-B10-329156	7
Fagron	Erythromycin	16G06-B03-324466	6
Fagron	Erythromycin	14/30-B01-303884	1
Fagron	Erythromycin	16J25-B10-329155	1
Fagron	Erythromycin	16J24-B10-329155	3
Fagron	Erythromycin	1430-B01-303884	1
Fagron	Erythromycin	15G07-B05-312356	2
Fagron	Erythromycin	1430B01303884	1
Fagron	Erythromycin	14B06-35-290873	1
Fagron	Erythromycin	14/30B0130050	1
Fagron	Erythromycin	13H22-N070126/246	1
Fagron	Erythromycin	13H22-N07	6
Fagron	Erythromycin	EB078B13	1
Fagron	Erythromycin	13A02-N-06	1
Fagron	Erythromycin	13D30-N01003/456	1
Fagron	Erythromycin	13F12-N21	25
Fagron	Erythromycin	12E24-No4	1
Fagron	Erythromycin	13D30N01	1
Fagron	Erythromycin	13D30-N010278/730	1
Fagron	Erythromycin	13D30-N010336/730	1
Fagron	Erythromycin	13d30-n01	1
Fagron	Erythromycin	13A02-N060571/680	1
Fagron	Erythromycin	13A02-No6	1
Fagron	Erythromycin	13D30-N01	26
Fagron	Erythromycin	13A02N06	1
Fagron	Erythromycin	E-030613-1	1
Fagron	Erythromycin	E-250213-2	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Erythromycin	11J06-No5	1
Fagron	Erythromycin	13A02-N06	11
Fagron	Erythromycin	12J30-N01	13
Fagron	Erythromycin	4212E-01362	1
Fagron	Erythromycin	510211362	1
Fagron	Erythromycin	11J06-N05	2
Fagron	Erythromycin	EBS/1201-008	2
Fagron	Erythromycin	12E24-N04	18
Fagron	Erythromycin	13H22-N070228/246	1
Fagron	Erythromycin	EB032J13	4
Fagron	Erythromycin	14B06-B35	11
Fagron	Erythromycin	14B06-B35-290875	5
Fagron	Erythromycin	14I30-B01-3000060	1
Fagron	Erythromycin	14I30-B01-300063	1
Fagron	Erythromycin	14F13-B12-300061	1
Fagron	Erythromycin	14I30B01	1
Fagron	Erythromycin	14/30B01300060	1
Fagron	Erythromycin	14i30b01-300059	1
Fagron	Erythromycin	14I30-B01-300060	2
Fagron	Erythromycin	14/30-B01-300060	3
Fagron	Erythromycin	14I30-B01-300062	2
Fagron	Erythromycin	EB022E14	2
Fagron	Erythromycin	14F13-B02-293511	2
Fagron	Erythromycin	14I30-B01	6
Fagron	Erythromycin	14I30-B01-303884	5
Fagron	Erythromycin	14F13-B02-293487	1
Fagron	Erythromycin	14B06-B35-293508	3
Fagron	Erythromycin	14F13-B02	3
Fagron	Erythromycin	14F13-B02-293490	2
Fagron	Erythromycin	14F13-B02-293488	3
Fagron	Erythromycin	14b06b5293508	1
Fagron	Erythromycin	14B06-B35-293509	6
Fagron	Erythromycin	14b06-b35-290872	1
Fagron	Erythromycin	13f12-n21	1
Fagron	Erythromycin	14B06-B35-290866	4
Fagron	Erythromycin	14B06-B35-290874	3
Fagron	Erythromycin	14B06-B35290873	1
Fagron	Erythromycin	14B06-B35-290873	3
Fagron	Erythromycin	13175810	3
Fagron	Erythromycin	16J24-B10-329158	3
Fagron	Erythromycin	14B06-B35-290872	5
Fagron	Erythromycin	18f01-b0a-353587	1
Fagron	Erythromycin	16J24-B10-339488	2
Fagron	Erythromycin	16J24-B10-3367363	1
Fagron	Erythromycin	16J24-B10-346206	3
Fagron	Erythromycin	16/24-B10-341875	1
Fagron	Erythromycin	18F05-B01-350907	1
Fagron	Erythromycin	18F05-B01-350906	1
Fagron	Erythromycin	CEPEB001E18	1
Fagron	Erythromycin	18F05-B01-350908	1
Fagron	Erythromycin	18F-05-B01	1
Fagron	Erythromycin	18F05-B01-352501	3
Fagron	Erythromycin	17088211	1
Fagron	Erythromycin	18F05-B01-353586	2
Fagron	Erythromycin	18f05-b01-350908	1
Fagron	Erythromycin	18f05-b01-352501	1
Fagron	Erythromycin	18F05-B01-352500	2
Fagron	Erythromycin	18F05-B01-353587	1
Fagron	Erythromycin	18f05-b01	1
Fagron	Erythromycin	700770-0002	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Erythromycin	16/24-B10-339487	2
Fagron	Erythromycin	16J24-B10-341875	5
Fagron	Erythromycin	16J24-b10-329155	1
Fagron	Erythromycin	16j24-b10	1
Fagron	Erythromycin	16j24-b10-329155	1
Fagron	Erythromycin	0EBA/0020516	2
Fagron	Erythromycin	16J24-B10-339487	5
Fagron	Escin soluble	30691/M3	1
Fagron	Estradiol benzoate	15C10-B04-305523	7
Fagron	Estradiol benzoate	13I11-N01	1
Fagron	Estradiol benzoate	15c10-b04-305523	1
Fagron	Estradiol benzoate	14B20-B02-294949	3
Fagron	Estradiol benzoate	17L01B05343714	1
Fagron	Estradiol benzoate	12H02-N02	3
Fagron	Estradiol benzoate	15E08-B08-311705	10
Fagron	Estradiol benzoate	12H02-N052	1
Fagron	Estradiol benzoate	81403680	1
Fagron	Estradiol benzoate	14B20-B02-2920001	1
Fagron	Estradiol benzoate	13/11-N01	2
Fagron	Estradiol benzoate	13I11N01	1
Fagron	Estradiol benzoate	81403690	2
Fagron	Estradiol benzoate	17L01-B05-3434714	1
Fagron	Estradiol benzoate	15E08-B08	1
Fagron	Estradiol benzoate	14B20-B02-292001	8
Fagron	Estradiol benzoate	12L14-N02	1
Fagron	Estradiol benzoate	14B20-B02	4
Fagron	Estradiol benzoate	13C15-B07-288645	1
Fagron	Estradiol benzoate	14b20-b02-294949	1
Fagron	Estradiol benzoate	11H01-N12	1
Fagron	Estradiol benzoate	12B08-N05	2
Fagron	Estradiol benzoate	17L01-B05-343714	1
Fagron	Estradiol benzoate	2x12L14-N02	1
Fagron	Estradiol benzoate	12D20N06	1
Fagron	Estradiol hemihydrate	15D01-B02-312068	2
Fagron	Estradiol hemihydrate	16h04-b01-334200	1
Fagron	Estradiol hemihydrate	16D12B01326783	1
Fagron	Estradiol hemihydrate	16H04-B01-328738	1
Fagron	Estradiol hemihydrate	16h04-b01-328738	1
Fagron	Estradiol hemihydrate	15D01-B02-319032	3
Fagron	Estradiol hemihydrate	18F18-B01-351790	1
Fagron	Estradiol hemihydrate	18F18-B01-351946	1
Fagron	Estradiol hemihydrate	16D12-B01-326783	5
Fagron	Estradiol hemihydrate	18F18-B01-351789	1
Fagron	Estradiol hemihydrate	9J27-N05	1
Fagron	Estradiol hemihydrate	15D01-B02-308568	1
Fagron	Estradiol hemihydrate	84201170	1
Fagron	Estradiol hemihydrate	14G24-B04-297404	1
Fagron	Estradiol hemihydrate	11A21-N01	1
Fagron	Estradiol hemihydrate	14G24-B04-297406	1
Fagron	Estradiol hemihydrate	14G24-B04-297405	1
Fagron	Estradiol hemihydrate	14G24-B04-297402	2
Fagron	Estradiol hemihydrate	13k22n01	1
Fagron	Estradiol hemihydrate	83100820	2
Fagron	Estradiol hemihydrate	13K22N01	1
Fagron	Estradiol hemihydrate	11031401	1
Fagron	Estradiol hemihydrate	13K22-N01	7
Fagron	Estradiol hemihydrate	12L29-N07	1
Fagron	Estradiol hemihydrate	13E07-N02	3
Fagron	Estradiol hemihydrate	9G07N01	1
Fagron	Estradiol hemihydrate	302446	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Estradiol hemihydrate	11H23-N02	1
Fagron	Estradiol hemihydrate	12L19-N07	6
Fagron	Estradiol hemihydrate	15D01-B02-308190	1
Fagron	Estradiol hemihydrate	12F29-N03	7
Fagron	Estradiol hemihydrate	1802013-01	5
Fagron	Estradiol hemihydrate	14G24-B04-302446	3
Fagron	Estradiol hemihydrate	16H04-B01-334200	1
Fagron	Estradiol hemihydrate	16H04-B01-339490	1
Fagron	Estradiol valerate	13g09-n05	1
Fagron	Estradiol valerate	16F09-B01-328135	3
Fagron	Estradiol valerate	84404140	3
Fagron	Estradiol valerate	15C10B05-326087	1
Fagron	Estradiol valerate	13I13-N02	1
Fagron	Estradiol valerate	15C10B05321675	1
Fagron	Estradiol valerate	14E26-B02-294845	5
Fagron	Estradiol valerate	14E26-B02-294847	1
Fagron	Estradiol valerate	15C10-B05-321675	5
Fagron	Estradiol valerate	19c10-b05-318405	1
Fagron	Estradiol valerate	14E26-B02-296348	6
Fagron	Estradiol valerate	14E26-B02-302873	3
Fagron	Estradiol valerate	14E26B02	2
Fagron	Estradiol valerate	15C10-B05-318405	2
Fagron	Estradiol valerate	15C10-B05	2
Fagron	Estradiol valerate	15C10-B05-314840	4
Fagron	Estradiol valerate	15c10-b05-309186	1
Fagron	Estradiol valerate	15C10-B05-305396	1
Fagron	Estradiol valerate	15C10-B05-309186	3
Fagron	Estradiol valerate	15C10-B05-308851	5
Fagron	Estradiol valerate	24052014a	1
Fagron	Estradiol valerate	16F09Bo1328135	1
Fagron	Estradiol valerate	13G09N05	1
Fagron	Estradiol valerate	15C10-B05_314840	1
Fagron	Estradiol valerate	11J13-No9	1
Fagron	Estradiol valerate	11K04-N04	1
Fagron	Estradiol valerate	12D18-N02	3
Fagron	Estradiol valerate	18G05-B02-353985	2
Fagron	Estradiol valerate	18g05-b02	1
Fagron	Estradiol valerate	12K05-N02	2
Fagron	Estradiol valerate	11k04-n04	1
Fagron	Estradiol valerate	15C10-B05-313568	1
Fagron	Estradiol valerate	13A07-N11	5
Fagron	Estradiol valerate	12G05N09	1
Fagron	Estradiol valerate	16F09-B01-332756	3
Fagron	Estradiol valerate	82201620	2
Fagron	Estradiol valerate	11G18-N04	1
Fagron	Estradiol valerate	13G09-N05	6
Fagron	Estradiol valerate	13D16N10	1
Fagron	Estradiol valerate	11J13-N09	1
Fagron	Estradiol valerate	16F09-B01_332756	1
Fagron	Estradiol valerate	13D16-N10	4
Fagron	Estradiol valerate	15C10-B05-305395	1
Fagron	Estriol	12F08-N06	1
Fagron	Estriol	13D19-N02	1
Fagron	Estriol	13H28-N01	3
Fagron	Estriol	12H06-No2	2
Fagron	Estriol	12H06-N02	6
Fagron	Estriol	81100490	1
Fagron	Estriol	13C15-N04	3
Fagron	Estriol	12K07-N15	2
Fagron	Estriol	86200800	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Estriol	16/19-B01-351346	1
Fagron	Estriol	16/19-B01-351346	1
Fagron	Estriol	16/19-B01-340536	1
Fagron	Estriol	16I19-B01-331467	1
Fagron	Estriol	16I9-B01	1
Fagron	Estriol	16/19-B01-332862	1
Fagron	Estriol	12A26-N02	2
Fagron	Estriol	16I19-B01-331468	1
Fagron	Estriol	16I19-B01-332862	1
Fagron	Estriol	16I19-B01	1
Fagron	Estriol	16D15-B03-323017	1
Fagron	Estriol	14E13-B01-294034	1
Fagron	Estriol	11D13-N11	1
Fagron	Estriol	15A15-B07-312064	1
Fagron	Estriol	14E13-B01-294033	1
Fagron	Estriol	13C15-N04	1
Fagron	Estriol	14b06b11291565	1
Fagron	Estriol	14E13-B01-294034	2
Fagron	Estriol	14B06-B11-294565	1
Fagron	Estriol	14B06-B11-291565	1
Fagron	Estriol	14B06-B11-291563	2
Fagron	Estriol	11J13-N07	3
Fagron	Estriol	16D15-B03-323016	2
Fagron	Ethacridine lactate m...	16E20-B11-328430	1
Fagron	Ethacridine lactate m...	12J24-N04	1
Fagron	Ethacridine lactate m...	11/07-N05	1
Fagron	Ethacridine lactate m...	12A02-N10	1
Fagron	Ethacridine lactate m...	12G02-N135-15	1
Fagron	Ethacridine lactate m...	12G02-N13	6
Fagron	Ethacridine lactate m...	201205001	2
Fagron	Ethacridine lactate m...	12C30-N07	1
Fagron	Ethacridine lactate m...	13/24-N03	4
Fagron	Ethacridine lactate m...	13124-N03	1
Fagron	Ethacridine lactate m...	12/24-N04	1
Fagron	Ethacridine lactate m...	14A15-B04-300493	6
Fagron	Ethacridine lactate m...	14A15-B04-300493	1
Fagron	Ethacridine lactate m...	16e20-b11-322514	1
Fagron	Ethacridine lactate m...	13/24-N03	1
Fagron	Ethacridine lactate m...	14A15-B04	1
Fagron	Ethacridine lactate m...	16E20-B11-322514	5
Fagron	Ethacridine lactate m...	12324N03	1
Fagron	Ethacridine lactate m...	12/24-N04	1
Fagron	Ethacridine lactate m...	13/24-N03	2
Fagron	Ethinylestradiol	16B17-B01-333326	1
Fagron	Ethinylestradiol	13G03-N04	1
Fagron	Ethinylestradiol	13I09-N04	1
Fagron	Ethinylestradiol	12D19-N04	1
Fagron	Ethinylestradiol	14e15b03296049	1
Fagron	Ethinylestradiol	84100720	1
Fagron	Ethinylestradiol	16B17-B01	1
Fagron	Ethinylestradiol	16B17-B01-320666	1
Fagron	Ethinylestradiol	18E25-B01-351832	2
Fagron	Ethinylestradiol	18E25-B01-353920	1
Fagron	Ethyl hydrogen fumara...	14H27-B04-299502	6
Fagron	Fampridine	18121801	1
Fagron	Fampridine	16/19-Bo4-330727	1
Fagron	Fampridine	13C07-B02-289395	1
Fagron	Fampridine	14D11-B06-296110	1
Fagron	Fampridine	14D11-B40-301641	1
Fagron	Fampridine	17/19-B01-339467	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Fampridine	16D11-B07	1
Fagron	Fampridine	17D19-B01-333745	1
Fagron	Fampridine	14D11-B40-304961	1
Fagron	Fenugreek, powdered	12149603	1
Fagron	Finasteride (form III)	19041707	1
Fagron	Finasteride (form III)	14111305	1
Fagron	Finasteride (form III)	17011303	1
Fagron	Finasteride (form III)	13I03-N01	2
Fagron	Finasteride (form III)	14L23-B02-306981	1
Fagron	Fluocinolonacetoneid, ...	13E22-N10	1
Fagron	Fluocinolonacetoneid, ...	13L30-B35-295489	1
Fagron	Fluocinolonacetoneid, ...	14L19-B04-615699	1
Fagron	Fluocinolonacetoneid, ...	14L19-B04	4
Fagron	Fluocinolonacetoneid, ...	13L30B35295489	1
Fagron	Fluocinolonacetoneid, ...	14L19-B04-317248	1
Fagron	Fluocinolonacetoneid, ...	15K06-B02	1
Fagron	Fluocinolonacetoneid, ...	15K06-B02-336730	2
Fagron	Folic acid hydrate	16C03-B06	1
Fagron	Folic acid hydrate	13A11-N01	1
Fagron	Folic acid hydrate	UT11110182	1
Fagron	Folic acid hydrate	13G16-N02	4
Fagron	Frankincense		1
Fagron	Fructose	16H05-B01	2
Fagron	Fructose	1521-B01-323021	1
Fagron	Fructose	14K17-B99	2
Fagron	Fructose	12H01-N03	2
Fagron	Fructose	13A28-N01	11
Fagron	Fructose	15/21-B01-312700	1
Fagron	Fructose	1312A-03394	1
Fagron	Fructose	15121-B01-312701	1
Fagron	Fructose	15A16-B07-308950	3
Fagron	Fructose	12Ho1-N03	1
Fagron	Fructose	15A16-B07-303124	1
Fagron	Fructose	12K28-N06	3
Fagron	Fructose	1170	1
Fagron	Fructose	13098906	5
Fagron	Fructose	15A16-B07	3
Fagron	Fructose	13E27-N07	3
Fagron	Fructose	13C25-N03	3
Fagron	Fructose	7K21-N12	1
Fagron	Fructose	13J18-B06	2
Fagron	Fructose	1107-N06	1
Fagron	Fructose	13J18-B06-286899	1
Fagron	Fructose	13J18-B06-286894	1
Fagron	Fructose	14E01-B04-293505	8
Fagron	Fructose	14E1-B04-293506	1
Fagron	Fructose	14E01-B04-29350	1
Fagron	Fructose	15A16-B07-305241	2
Fagron	Fructose	13E27-n07	1
Fagron	Fructose	14E01-B04	3
Fagron	Fructose	14E01-B04-297098	2
Fagron	Fructose	14k17-B99	2
Fagron	Fructose	14E01B04293506	1
Fagron	Fructose	15121-B01-323021	1
Fagron	Fructose	508Q-03394	1
Fagron	Fructose	FR1443072	8
Fagron	Fructose	15I21-B01-323020	1
Fagron	Fructose	15/21-B01	1
Fagron	Fructose	16H05-B01-326259	1
Fagron	Fructose	3014Q-03394	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Fructose	15/21-B01-323021	1
Fagron	Fructose	15I21-B01-323021	5
Fagron	Fructose	326259	1
Fagron	Fructose	16H05-B01-330164	2
Fagron	Fructose	12A26-N04	1
Fagron	Fructose	516Q-03394	1
Fagron	Fructose	12E29-N08	3
Fagron	Fructose	15I21-B01	2
Fagron	Fuchsine	12A03-N06	1
Fagron	Fuchsine	12F28-N02	1
Fagron	Fuchsine	12G18-N01	1
Fagron	Fuchsine	13A30-N05	1
Fagron	Fuchsine	13H23-M04	2
Fagron	Furosemide	13E28-B40-289433	1
Fagron	Furosemide	13E28-B40-295538	1
Fagron	Furosemide	13E28-B40	2
Fagron	Furosemide	12I13-No6	1
Fagron	Furosemide	1095H2RI	1
Fagron	Furosemide	13A22-N08	3
Fagron	Furosemide	11K21-N26	1
Fagron	Furosemide	13A22-N08	1
Fagron	Furosemide	15L10-B02-316506	2
Fagron	Furosemide	17d06-b01-335672	1
Fagron	Furosemide	AM20216066	1
Fagron	Furosemide	16/26-B04-329253	1
Fagron	Furosemide	17D06-B01-335672	1
Fagron	Furosemide	18C27-B02	1
Fagron	Furosemide	12A18-N03	1
Fagron	Furosemide	14K10-B02-301746	1
Fagron	Fusidic acid	FS655	1
Fagron	Fusidic acid	13A07-N10	2
Fagron	Fusidic acid	13H19-N05	1
Fagron	Fusidic acid	12G04-No1	1
Fagron	Fusidic acid	13L13-B01-289514	1
Fagron	Fusidic acid	17G06-B99-337953	1
Fagron	Fusidic acid	12B23-N03	1
Fagron	Fusidic acid	13L13-B01-289513	2
Fagron	Fusidic acid	4113U-07034	1
Fagron	Fusidic acid	14E01-B08-299791	3
Fagron	Fusidic acid	15/10-B02-314192	1
Fagron	Fusidic acid	13L13-B01	1
Fagron	Fusidic acid	1510-B02_319231	1
Fagron	Fusidic acid	15I10B02-319231	1
Fagron	Fusidic acid	15I10-B02-326182	1
Fagron	Fusidic acid	15I10-B02	1
Fagron	Fusidic acid	17J06-B02-343677	1
Fagron	Fusidic acid	116E-03682	1
Fagron	Fusidic acid	15/10-b02	2
Fagron	Gentamicin sulfate	11/09-N10	1
Fagron	Gentamicin sulfate	12D26-N05	8
Fagron	Gentamicin sulfate	12D6-N05	1
Fagron	Gentamicin sulfate	16829-B01-328731	1
Fagron	Gentamicin sulfate	13D24-N05	7
Fagron	Gentamicin sulfate	FG1212646	1
Fagron	Gentamicin sulfate	13D24-N05	1
Fagron	Gentamicin sulfate	13I19-B05	1
Fagron	Gentamicin sulfate	14A30-B04-290861	3
Fagron	Gentamicin sulfate	13I19-B05-286568	1
Fagron	Gentamicin sulfate	FG1305122	1
Fagron	Gentamicin sulfate	18G06-B09-354797	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Gentamicin sulfate	13/19-B05-286568	1
Fagron	Gentamicin sulfate	10A29-N06	1
Fagron	Gentamicin sulfate	5214A-01582	1
Fagron	Gentamicin sulfate	11K21-N17	2
Fagron	Gentamicin sulfate	11k21-n17	1
Fagron	Gentamicin sulfate	18G06-B09-351916	1
Fagron	Gentamicin sulfate	18G06-B09	1
Fagron	Gentamicin sulfate	18C12B02350204	1
Fagron	Gentamicin sulfate	18C12-B02-350204	2
Fagron	Gentamicin sulfate	18C12-B02	2
Fagron	Gentamicin sulfate	18C12-B02-350	1
Fagron	Gentamicin sulfate	12117-no2	1
Fagron	Gentamicin sulfate	18C12-B02-349025	5
Fagron	Gentamicin sulfate	15121722	2
Fagron	Gentamicin sulfate	17G27-B05	1
Fagron	Gentamicin sulfate	16B29-B01-328915	1
Fagron	Gentamicin sulfate	12B21-N02	2
Fagron	Gentamicin sulfate	13/19-B05	1
Fagron	Gentamicin sulfate	1221-n02	1
Fagron	Gentamicin sulfate	17G27-B05-340121	1
Fagron	Gentamicin sulfate	14A30-B04	1
Fagron	Gentamicin sulfate	1319-B05-286565	1
Fagron	Gentamicin sulfate	14A30-B04	2
Fagron	Gentamicin sulfate	14K28-B06-301951	1
Fagron	Gentamicin sulfate	14K28-B06	4
Fagron	Gentamicin sulfate	14K28-B06-303640	3
Fagron	Gentamicin sulfate	14A30-B04-290858	1
Fagron	Gentamicin sulfate	15G29-B09-326635	1
Fagron	Gentamicin sulfate	15G29B09322936	1
Fagron	Gentamicin sulfate	10C22N02	1
Fagron	Gentamicin sulfate	15G29-B09-322936	2
Fagron	Gentamicin sulfate	15G29-B09-326632	4
Fagron	Gentamicin sulfate	15G29-B09-322945	1
Fagron	Gentamicin sulfate	15E04-B05-319740	1
Fagron	Gentamicin sulfate	15803-B07-309482	1
Fagron	Gentamicin sulfate	11604054-01	1
Fagron	Gentamicin sulfate	15G29-B09-320305	3
Fagron	Gentamicin sulfate	15E04-B05-316255	2
Fagron	Gentamicin sulfate	15B03-B07	1
Fagron	Gentamicin sulfate	15B03-B07-309482	2
Fagron	Gentamicin sulfate	14K28B06303646	2
Fagron	Gentamicin sulfate	15b03-b07-309482	1
Fagron	Gentamicin sulfate	14k28-b06-303640	2
Fagron	Gentamicin sulfate	14k28-B06-303640	1
Fagron	Gentamicin sulfate	12D26-N05	2
Fagron	Gentamicin sulfate	14e08b03300619	1
Fagron	Gentamicin sulfate	14081307	3
Fagron	Gentamicin sulfate	15G29-B09-322928	1
Fagron	Gentamicin sulfate	16B29-B01-328731	4
Fagron	Gentamicin sulfate	14E08-B03-294848	3
Fagron	Gentamicin sulfate	15G29B0932632	1
Fagron	Gentamicin sulfate	13F13-N01	1
Fagron	Gentamicin sulfate	14a30-b04-290858	1
Fagron	Gentamicin sulfate	13D24-N05	1
Fagron	Gentamicin sulfate	14E08-B03	2
Fagron	Gentamicin sulfate	14E08-B03-294835	2
Fagron	Gentamicin sulfate	19-B05-286568	1
Fagron	Gentamicin sulfate	14E08-B03-300619	2
Fagron	Glucose	14142303	3
Fagron	Glucose monohydrate	13H14-B06-289483	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Glucose monohydrate	15B01291137	1
Fagron	Glucose monohydrate	1521-B06-309508	1
Fagron	Glucose monohydrate	14B28-B02	2
Fagron	Glucose monohydrate	14I04B8298602	1
Fagron	Glucose monohydrate		2
Fagron	Glucose monohydrate	E0AJ5	1
Fagron	Glucose monohydrate	15A21-B06-309508	3
Fagron	Glucose monohydrate	1407-808-307707	1
Fagron	Glucose monohydrate	14B04-B01-291138	2
Fagron	Glucose monohydrate	14D30-B06-294678	1
Fagron	Glucose monohydrate	13C25N09	2
Fagron	Glucose monohydrate	14B28-B02-291354	3
Fagron	Glucose monohydrate	13B01-N03	5
Fagron	Glucose monohydrate	13A25-n01	1
Fagron	Glucose monohydrate	14I04-B08-298601	1
Fagron	Glucose monohydrate	4562135	1
Fagron	Glucose monohydrate	3413E-02209	1
Fagron	Glucose monohydrate	1404-B08-298601	2
Fagron	Glucose monohydrate	4213Q-02209	1
Fagron	Glucose monohydrate	1404-B08-298604	1
Fagron	Glucose monohydrate	14104-B08-302707	1
Fagron	Glucose monohydrate	14I04-808	1
Fagron	Glucose monohydrate	14B04-B01-291137	8
Fagron	Glucose monohydrate	14B28-B02-291353	10
Fagron	Glucose monohydrate	14B04-B01-291143	2
Fagron	Glucose monohydrate	14I04-B08-298602	3
Fagron	Glucose monohydrate	14B04-B01	4
Fagron	Glucose monohydrate	8G31-N09	1
Fagron	Glucose monohydrate	151780-0004	1
Fagron	Glucose monohydrate	E571X	1
Fagron	Glucose monohydrate	13A25-N01	2
Fagron	Glucose monohydrate	13A25-No1	1
Fagron	Glucose monohydrate	13A_N01	1
Fagron	Glucose monohydrate	12B09-N05	1
Fagron	Glucose monohydrate	13A25-N01	12
Fagron	Glucose monohydrate	270313PML01	1
Fagron	Glucose monohydrate	12G30-N01	2
Fagron	Glucose monohydrate	12G30-N01	13
Fagron	Glucose monohydrate	E524W	1
Fagron	Glucose monohydrate	13C25-N09	9
Fagron	Glucose monohydrate	14I04-B08	5
Fagron	Glutamic acid	3031502	1
Fagron	Glycerol monostearate...	15280903	3
Fagron	Glycerol monostearate...	7G19-N09	1
Fagron	Glycerol monostearate...	12K05-N06	1
Fagron	Glycerol monostearate...	12A27-N11	1
Fagron	Glycopyrronium bromide	18C30-B07-350744	2
Fagron	Gramicidin	11F-24-N03	1
Fagron	Gramicidin	12K22-N12	1
Fagron	Gramicidin	17E09-B02-341211	1
Fagron	Gramicidin	14H06-B03-297614	2
Fagron	Gramicidin	B3200266	1
Fagron	Gramicidin	15C11-B04-306766	1
Fagron	Griseofulvin	18L04-B01-357660	1
Fagron	Griseofulvin	16F13-B06-323067	3
Fagron	Griseofulvin	17-L05-352461	1
Fagron	Griseofulvin	15B10-B06-306547	9
Fagron	Griseofulvin	15810-B06	1
Fagron	Griseofulvin	16f13-b06-326074	1
Fagron	Griseofulvin	15L03-B04-315624	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Griseofulvin	14E16-B01-300367	1
Fagron	Griseofulvin	15L03-b04-315624	1
Fagron	Griseofulvin	17L05-B05-352461	1
Fagron	Griseofulvin	18L04-B01-357660	1
Fagron	Griseofulvin	15L03-B315624	1
Fagron	Griseofulvin	15B10-B06-306663	1
Fagron	Griseofulvin	17L005-B05-352461	1
Fagron	Griseofulvin	11J11-N13	1
Fagron	Griseofulvin	12F18-N04	1
Fagron	Griseofulvin	16F13-B06-326073	2
Fagron	Griseofulvin	16F13-B06-330832	4
Fagron	Griseofulvin	16F13-B06-330833	4
Fagron	Griseofulvin	15B10-B06	1
Fagron	Griseofulvin	16F13B06330833	1
Fagron	Griseofulvin	15B10-B06-396547	1
Fagron	Griseofulvin	16f13-b06-330833	1
Fagron	Griseofulvin	16F13-B06	1
Fagron	Griseofulvin	HH160510	1
Fagron	Griseofulvin	17L05-B05-34999	1
Fagron	Griseofulvin	17L05-B05	1
Fagron	Griseofulvin	HH171022	1
Fagron	Griseofulvin	16F13-B06-323068	3
Fagron	Griseofulvin	15062015A	1
Fagron	Griseofulvin	14E16-B01-300366	3
Fagron	Griseofulvin	17L05B05352461	2
Fagron	Griseofulvin	13L13-B08	1
Fagron	Griseofulvin	10F03-N12	1
Fagron	Griseofulvin	13L13-B08-290316	3
Fagron	Griseofulvin	12D04-N12	1
Fagron	Griseofulvin	130922	2
Fagron	Griseofulvin	13F12-N22	1
Fagron	Griseofulvin	17L05-B05-351232	1
Fagron	Griseofulvin	13D19-N03	1
Fagron	Griseofulvin	13L-13B08	1
Fagron	Hard fat in pastilles	16H01-B03-334976	1
Fagron	Hard fat in pastilles	14L01-B06-303380	3
Fagron	Hard fat in pastilles	15F03-B06-308708	3
Fagron	Hard fat in pastilles	15F03-B06	4
Fagron	Hard fat in pastilles	14L01-B06	1
Fagron	Hard fat in pastilles	12G04-N18	1
Fagron	Hard fat in pastilles	12I04-N11	1
Fagron	Hard fat in pastilles	606469	1
Fagron	Hard fat in pastilles	13D10-N06	1
Fagron	Hard fat in pastilles	14G03-B03-296731	1
Fagron	Hard fat W35	17k09-b03-351980	1
Fagron	Hard fat W35	17K09-B03-351980	2
Fagron	Hard fat W35	17K09-B03-346709	1
Fagron	Hard fat W35	15J08-B04-322726	1
Fagron	Hard fat W35	17K09-B03-348770	3
Fagron	Hard fat W35	16K08-B10-339077	3
Fagron	Hard fat W35	3914E-02328	1
Fagron	Hard fat W35	15J08-B04-322727	1
Fagron	Hard fat W35	15J08-B04-324217	1
Fagron	Hard fat W35	15J08-B04-324959	2
Fagron	Hard fat W35	16K08-B10-332763	2
Fagron	Hard fat W35	15A07-B12-309245	1
Fagron	Heavy alkaline bismut...	13B25-N02	3
Fagron	Heavy alkaline bismut...	1001311186	1
Fagron	Heavy alkaline bismut...	11/13-N01	1
Fagron	Heavy alkaline bismut...	11D28-N16	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Heavy alkaline bismut...	13C15-N05	3
Fagron	Heavy alkaline bismut...	12C08-N06	3
Fagron	Heavy alkaline bismut...	13H28-B04-291385	1
Fagron	Heavy alkaline bismut...	1309211186	1
Fagron	Heavy alkaline bismut...	12D26-N03	1
Fagron	Heparin sodium	14D09-B08-293944	1
Fagron	Heparin sodium	12G19-B02-313422	1
Fagron	Heparin sodium	14D09-B08-307077	1
Fagron	Heparin sodium	13G03-N13	3
Fagron	Heparin sodium	11190188SH	2
Fagron	Heparin sodium	12H30-N07	2
Fagron	Heparin sodium	12120-B01-323107	1
Fagron	Heparin sodium	12190453CQ	1
Fagron	Heparin sodium	16I26-B02	1
Fagron	Heparin sodium	15A20-B02-307307	1
Fagron	Heparin sodium	16I26-B02-332108	1
Fagron	Heparin sodium	16I26-B02-327885	2
Fagron	Histamine dihydrochlo...	L17100130-0F-222030	1
Fagron	Histamine dihydrochlo...	15E08-B06-313037	3
Fagron	Histamine dihydrochlo...	14J06B01299179	1
Fagron	Histamine dihydrochlo...	12E29-N09	3
Fagron	Histamine dihydrochlo...	12D02-N16	1
Fagron	Histamine dihydrochlo...	454545	1
Fagron	Histamine dihydrochlo...	12K06-N08	7
Fagron	Histamine dihydrochlo...	28051305	1
Fagron	Histamine dihydrochlo...	11E25-N07	2
Fagron	Hyaluronic acid sodiu...	L1311012-01	1
Fagron	Hydrochlorothiazide	12A03-B40	3
Fagron	Hydrochlorothiazide	14B10-B01-300822	2
Fagron	Hydrochlorothiazide	11L01-N02	2
Fagron	Hydrochlorothiazide	13008504	1
Fagron	Hydrochlorothiazide	11Lo1-No2	1
Fagron	Hydrochlorothiazide	14b10-b01-294558	1
Fagron	Hydrochlorothiazide	11L01-No2	1
Fagron	Hydrochlorothiazide	12L14-N11	4
Fagron	Hydrochlorothiazide	13C25-N08	6
Fagron	Hydrochlorothiazide	14B10-B01-294558	6
Fagron	Hydrochlorothiazide	11E04-N10	1
Fagron	Hydrochlorothiazide	2002311234	1
Fagron	Hydrochlorothiazide	12A03-B40-286110	1
Fagron	Hydrochlorothiazide	CQ1108	1
Fagron	Hydrochlorothiazide	12A03-B40-286111	3
Fagron	Hydrochlorothiazide	CR0601	1
Fagron	Hydrochlorothiazide	5013Q-03660	1
Fagron	Hydrochlorothiazide	11G12-N12	1
Fagron	Hydrochlorothiazide	15K16-B01-315836	2
Fagron	Hydrochlorothiazide	14B10B01	1
Fagron	Hydrochlorothiazide	17L20-B02-355636	2
Fagron	Hydrochlorothiazide	14B10-B01-305043	4
Fagron	Hydrochlorothiazide	17L20-B02-349593	1
Fagron	Hydrochlorothiazide	17L20B05	1
Fagron	Hydrochlorothiazide	17L20-B02-347789	3
Fagron	Hydrochlorothiazide	17L20-B02-345079	3
Fagron	Hydrochlorothiazide	17H28-B01-340561	1
Fagron	Hydrochlorothiazide	14B10-B01	2
Fagron	Hydrochlorothiazide	16F17-B07-324062	6
Fagron	Hydrochlorothiazide	16F17-B07	2
Fagron	Hydrochlorothiazide	11101-n02	2
Fagron	Hydrochlorothiazide	Co1-20150412	1
Fagron	Hydrochlorothiazide	15K16-B01	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hydrochlorothiazide	15E04-B03-309670	4
Fagron	Hydrochlorothiazide	15E04-B03	1
Fagron	Hydrochlorothiazide	C01-20150412	1
Fagron	Hydrochlorothiazide	15K16-B01-319833	8
Fagron	Hydrochlorothiazide	15K16	1
Fagron	Hydrocortisone	14B04-B02-293423	3
Fagron	Hydrocortisone	13H08-B01-285473	1
Fagron	Hydrocortisone	13K22-B01-288284	4
Fagron	Hydrocortisone	13H08B01285473	1
Fagron	Hydrocortisone	K05B20131101	1
Fagron	Hydrocortisone	13K22-B01-288283	1
Fagron	Hydrocortisone	14B04-B02	2
Fagron	Hydrocortisone	19051403	1
Fagron	Hydrocortisone	14B04-B02-293424	4
Fagron	Hydrocortisone	14105-B02-300228	1
Fagron	Hydrocortisone	5061304	1
Fagron	Hydrocortisone	14B04-B02-293422	2
Fagron	Hydrocortisone	14/05-B02-300227	2
Fagron	Hydrocortisone	12I25-N01	2
Fagron	Hydrocortisone	12/25-N01	1
Fagron	Hydrocortisone	12G19-N01	2
Fagron	Hydrocortisone	11K18-N01	2
Fagron	Hydrocortisone	120601801	1
Fagron	Hydrocortisone	PH132	1
Fagron	Hydrocortisone	13H08-B01-285472	1
Fagron	Hydrocortisone	13H08-B01-28571	1
Fagron	Hydrocortisone	13K22-B01-288282	1
Fagron	Hydrocortisone	17E-10-B04-334857	1
Fagron	Hydrocortisone	17I12-B05	1
Fagron	Hydrocortisone	15/25-B01-318727	1
Fagron	Hydrocortisone	15125-B01-318486	1
Fagron	Hydrocortisone	17I12-B05-344220	1
Fagron	Hydrocortisone	17/12-B05-344220	2
Fagron	Hydrocortisone	16H01-B01	1
Fagron	Hydrocortisone	17E10B04335656	1
Fagron	Hydrocortisone	17G25-B02-340864	1
Fagron	Hydrocortisone	16A19-B06	1
Fagron	Hydrocortisone	17E10-B04-334858	3
Fagron	Hydrocortisone	14105-B02-300227	1
Fagron	Hydrocortisone	16j26-b04-330339	1
Fagron	Hydrocortisone	16H01-B01-324944	5
Fagron	Hydrocortisone	16A19-B06-321982	1
Fagron	Hydrocortisone	17E10-B04-336436	2
Fagron	Hydrocortisone	17E10-B04-336435	2
Fagron	Hydrocortisone	16A19-B06-320970	20
Fagron	Hydrocortisone	K05M20170702	1
Fagron	Hydrocortisone	17-12B05-342611	1
Fagron	Hydrocortisone	16A19-B06-321981	2
Fagron	Hydrocortisone	15I04-B03-315021	1
Fagron	Hydrocortisone	17E10-B04-334857	5
Fagron	Hydrocortisone	14 05-B02-300227	1
Fagron	Hydrocortisone	14I05-B02-300229	1
Fagron	Hydrocortisone	14I05-B02	2
Fagron	Hydrocortisone	17I12-B05-34554	1
Fagron	Hydrocortisone	16J26-B04-330340	1
Fagron	Hydrocortisone	17E10-B04	1
Fagron	Hydrocortisone	17/12-B05-344220	1
Fagron	Hydrocortisone	16H01-B01-324891	1
Fagron	Hydrocortisone	14/05-B02-300227	1
Fagron	Hydrocortisone	17C13-B06-333045	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hydrocortisone	15D14-B02	1
Fagron	Hydrocortisone	15/25-B01-313084	1
Fagron	Hydrocortisone	15D14-B02-311696	5
Fagron	Hydrocortisone acetate	16H18-B04-325093	2
Fagron	Hydrocortisone acetate	16J28-B05-330083	1
Fagron	Hydrocortisone acetate	11/01-N05	1
Fagron	Hydrocortisone acetate	12K15-N04	6
Fagron	Hydrocortisone acetate	16H18-B04-327325	1
Fagron	Hydrocortisone acetate	1112031	1
Fagron	Hydrocortisone acetate	16H18-B04	1
Fagron	Hydrocortisone acetate	13G09-N08	7
Fagron	Hydrocortisone acetate	7C10-B01-333961	1
Fagron	Hydrocortisone acetate	17L19	2
Fagron	Hydrocortisone acetate	17L19B04345374	1
Fagron	Hydrocortisone acetate	17G03-B06-345375	3
Fagron	Hydrocortisone acetate	17L19-B04-350448	1
Fagron	Hydrocortisone acetate	17L19-B04-348983	1
Fagron	Hydrocortisone acetate	17G03-B06-341304	1
Fagron	Hydrocortisone acetate	17L19-B04	1
Fagron	Hydrocortisone acetate	17086211	6
Fagron	Hydrocortisone acetate	18H07-B03-352228	2
Fagron	Hydrocortisone acetate	17G03-B06-341843	1
Fagron	Hydrocortisone acetate	17C10-B01-335774	1
Fagron	Hydrocortisone acetate	16C15-B04-321459	2
Fagron	Hydrocortisone acetate	17C10	2
Fagron	Hydrocortisone acetate	17C10-B01-333961	1
Fagron	Hydrocortisone acetate	17G03-B06-338750	1
Fagron	Hydrocortisone acetate	17C10-B01-333962	3
Fagron	Hydrocortisone acetate	18H07-B03-356259	1
Fagron	Hydrocortisone acetate	NCEHAC151101	1
Fagron	Hydrocortisone acetate	16J28-B05	1
Fagron	Hydrocortisone acetate	13J31-B07-290062	2
Fagron	Hydrocortisone acetate	15A03-B09-307133	1
Fagron	Hydrocortisone acetate	13J31-B07-290061	2
Fagron	Hydrocortisone acetate	HC0024	1
Fagron	Hydrocortisone acetate	15A30-B09-307133	4
Fagron	Hydrocortisone acetate	15A03-B09-307300	1
Fagron	Hydrocortisone acetate	13F25-B05-300290	2
Fagron	Hydrocortisone acetate	15H24-B01-314236	3
Fagron	Hydrocortisone acetate	13E31-N02	2
Fagron	Hydrocortisone acetate	14043906	7
Fagron	Hydrocortisone acetate	15H24-B01-316480	1
Fagron	Hydrocortisone acetate	14A13-B04-294158	2
Fagron	Hydrocortisone acetate	13E31-No2	1
Fagron	Hydrocortisone acetate	M39171	1
Fagron	Hydrocortisone acetate	8563642	1
Fagron	Hydrocortisone acetate	15H24-B01	1
Fagron	Hydrocortisone acetate	14A13-B04	1
Fagron	Hydrocortisone acetate	13J31-B07	3
Fagron	Hydrocortisone acetate	3613A-01718	1
Fagron	Hydrocortisone acetate	15J28-B12-317771	1
Fagron	Hydrocortisone acetate	13J31B07290061	1
Fagron	Hydrocortisone acetate	15h24-b01-314236	1
Fagron	Hydrocortisone acetate	13155311	1
Fagron	Hydrocortisone acetate	13A03-N04	4
Fagron	Hydrocortisone acetate	13/31-B07-290061	1
Fagron	Hydrocortisone acetate	13B25-N04	2
Fagron	Hydrocortisone butyrate	15G08-B16-314883	1
Fagron	Hydrocortisone butyrate	16j21-b01	1
Fagron	Hydroquinone		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hydroxyethylcellulose...	15E07-B05-316511	2
Fagron	Hydroxyethylcellulose...	14/04-B05-300231	1
Fagron	Hydroxyethylcellulose...	15E07-B05-316511	1
Fagron	Hydroxyethylcellulose...	14F04-B06-297846	1
Fagron	Hydroxyethylcellulose...	14I04-B05-300231	1
Fagron	Hydroxyethylcellulose...	15B04-B03-312859	4
Fagron	Hydroxyethylcellulose...	14i04-b05-309049	1
Fagron	Hydroxyethylcellulose...	K1373	1
Fagron	Hydroxyethylcellulose...	11I02-N02	1
Fagron	Hydroxyethylcellulose...	12A09-N06	2
Fagron	Hydroxyethylcellulose...	14I04-B05-309049	1
Fagron	Hydroxyethylcellulose...	14I04-B05-300231	1
Fagron	Hydroxyethylcellulose...	1404-B05-300231	1
Fagron	Hydroxyethylcellulose...	14Io4-B05	1
Fagron	Hydroxyethylcellulose...	14I04-B05-300224	1
Fagron	Hydroxyethylcellulose...	12C27-N01	1
Fagron	Hydroxyethylcellulose...	15086407	1
Fagron	Hydroxyethylcellulose...	15086403	1
Fagron	Hydroxyethylcellulose...	15H28-B10	2
Fagron	Hydroxyethylcellulose...	15H28-B10-319028	2
Fagron	Hydroxyethylcellulose...	16107-b04-331244	1
Fagron	Hydroxyethylcellulose...	16H05B05325803	1
Fagron	Hydroxyethylcellulose...	18C05-B05-352498	2
Fagron	Hydroxyethylcellulose...	11I02-N12	1
Fagron	Hydroxyethylcellulose...	18C05-B05-354803	2
Fagron	Hydroxyethylcellulose...	19C04-B08-191454	1
Fagron	Hydroxyethylcellulose...	12H02-N06	3
Fagron	Hydroxyethylcellulose...	13D26-N11	5
Fagron	Hydroxyethylcellulose...	14D04-B02-293558	1
Fagron	Hydroxyethylcellulose...	14I04-B05-300224	1
Fagron	Hydroxyethylcellulose...	15B04-B03-315495	1
Fagron	Hydroxyethylcellulose...	14I04-B05	1
Fagron	Hydroxyethylcellulose...	13H22-N11	8
Fagron	Hydroxyethylcellulose...	13L04-B03-289115	3
Fagron	Hydroxyethylcellulose...	13L04-B03	2
Fagron	Hydroxyethylcellulose...	14F04-B06-295130	2
Fagron	Hydroxyethylcellulose...	14f04-b06-295130	1
Fagron	Hydroxyethylcellulose...	14D04B02293555	1
Fagron	Hydroxyethylcellulose...	K1373	1
Fagron	Hydroxyethylcellulose...	17J05-B03-347890	7
Fagron	Hydroxyethylcellulose...	16E09-B05-321320	4
Fagron	Hydroxyethylcellulose...	17J05-B03-343927	1
Fagron	Hydroxyethylcellulose...	16L07-B04-331846	3
Fagron	Hydroxyethylcellulose...	14D04-B02-293558	1
Fagron	Hydroxyethylcellulose...	16H05-B05	1
Fagron	Hydroxyethylcellulose...	16E09-B05-327042	1
Fagron	Hydroxyethylcellulose...	16L07-B04-331244	2
Fagron	Hydroxyethylcellulose...	16H05-B05-325803	1
Fagron	Hydroxyethylcellulose...	17G05-B04-341204	11
Fagron	Hydroxyethylcellulose...	14D04-B02	2
Fagron	Hydroxyethylcellulose...	17G05-B04	3
Fagron	Hydroxyethylcellulose...	J0696	2
Fagron	Hydroxyethylcellulose...	17J05-B03-342018	2
Fagron	Hydroxyethylcellulose...	17j05B03342018	1
Fagron	Hydroxyethylcellulose...	17j05-b03-343927	1
Fagron	Hydroxyethylcellulose...	14352504	1
Fagron	Hydroxyethylcellulose...	130902002	1
Fagron	Hydroxyethylcellulose...	13d26-n16	1
Fagron	Hydroxyethylcellulose...	11I02-N02	1
Fagron	Hyoscine butylbromide	16I20-B08-326942	12

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hyoscine hydrobromide	14J27-B05-307020	1
Fagron	Hyoscine hydrobromide	18B23-B04-357420	1
Fagron	Hyoscine hydrobromide	14F11-B02-296840	1
Fagron	Hypromellose 4000	18A17-B04-346208	1
Fagron	Hypromellose 4000	12A30-N02	3
Fagron	Hypromellose 4000	15G08-B14	2
Fagron	Hypromellose 4000	15G08-B14-314845	8
Fagron	Hypromellose 4000	15g08-b14-314845	1
Fagron	Hypromellose 4000	15G08B14314185	2
Fagron	Hypromellose 4000	3G08012N11	1
Fagron	Hypromellose 4000	15g08-b14-310845	2
Fagron	Hypromellose 4000	15G08-B14-310845	8
Fagron	Hypromellose 4000	14g02-b11-302118	2
Fagron	Hypromellose 4000	15B02-B05	1
Fagron	Hypromellose 4000	1J31012N11	3
Fagron	Hypromellose 4000	14G02-B11-302119	1
Fagron	Hypromellose 4000	13H09-N04	7
Fagron	Hypromellose 4000	13D03-N03	4
Fagron	Hypromellose 4000	14G02-B11-302121	1
Fagron	Hypromellose 4000	1K01012N11	1
Fagron	Hypromellose 4000	14G02-B11	4
Fagron	Hypromellose 4000	14G02-B11-302118	6
Fagron	Hypromellose 4000	13E02-B03-290532	1
Fagron	Hypromellose 4000	2H12012N12	1
Fagron	Hypromellose 4000	14C-05B03	1
Fagron	Hypromellose 4000	14C05B03	1
Fagron	Hypromellose 4000	14G02-B11-296851	1
Fagron	Hypromellose 4000	14c05-B03-294708	1
Fagron	Hypromellose 4000	14C05-B03-294711	1
Fagron	Hypromellose 4000	700182-0001	2
Fagron	Hypromellose 4000	14C05-B03-294708	5
Fagron	Hypromellose 4000	13E02-B03-290531	2
Fagron	Hypromellose 4000	D180F8Q011	2
Fagron	Hypromellose 4000	16B16-B07-318276	1
Fagron	Hypromellose 4000	16B16-B07-318531	8
Fagron	Hypromellose 4000	18D04-B03-35144	1
Fagron	Hypromellose 4000	20110612-2	1
Fagron	Hypromellose 4000	18D04-B03-352746	1
Fagron	Hypromellose 4000	18D04-B03-350745	1
Fagron	Hypromellose 4000	10G09-N03	1
Fagron	Hypromellose 4000	18D04-B03	2
Fagron	Hypromellose 4000	12E02-N03	3
Fagron	Hypromellose 4000	349970	1
Fagron	Hypromellose 4000	12G23-N03	1
Fagron	Hypromellose 4000	ZH19012N12	1
Fagron	Hypromellose 4000	12G23-N03	1
Fagron	Hypromellose 4000	325994	1
Fagron	Hypromellose 4000	13H09-N04	3
Fagron	Hypromellose 4000	16B16-B07-325280	2
Fagron	Hypromellose 4000	12H20-N05	1
Fagron	Hypromellose 4000	13D03-N03	3
Fagron	Hypromellose 4000	13D03N03	1
Fagron	Hypromellose 4000	16B16-B07-321327	5
Fagron	Hypromellose 4000	16B-16-B07-321327	1
Fagron	Hypromellose 4000	YI16012N11	1
Fagron	Hypromellose 4000	16H23-B02-332066	6
Fagron	Hypromellose 4000	16B16-B07-321328	1
Fagron	Hypromellose 4000	16B16-B07-321579	1
Fagron	Hypromellose 4000	16B16_B07	1
Fagron	Hypromellose 4000	16B16-B07-323114	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Hypromellose 4000	12D03-N03	1
Fagron	Hypromellose 4000	13h09-m04	1
Fagron	Hypromellose 4000	1A01012N11	1
Fagron	Hypromellose 4000	14C05-B03	1
Fagron	Indigo carmine	11F09-N06	1
Fagron	Indigo carmine	12K06-N11	1
Fagron	Indigo carmine	12A20-N067-16	1
Fagron	Indometacin		1
Fagron	Indometacin	14I08-B06-343286	1
Fagron	Indometacin	17I08B06	1
Fagron	Indometacin	17I08-B06-339738	1
Fagron	Inositol	L13090200	1
Fagron	Inositol hexanicotina...	1811046001	1
Fagron	Iodoform	12L11-B41-290055	1
Fagron	Iodoform	12E30-N01	2
Fagron	Iodoform	18091803	1
Fagron	Ipratropium bromide	14B05-B16-294566	1
Fagron	Isosorbide dinitrate ...	17C29-B01_350469	1
Fagron	Isosorbide dinitrate ...	17c29-b01-350469	1
Fagron	Isosorbide dinitrate ...	17C29-B01-339448	2
Fagron	Isosorbide dinitrate ...	17c29-b01	1
Fagron	Isosorbide dinitrate ...	17c29-b01-348859	2
Fagron	Isosorbide dinitrate ...	17C29-801-347885	1
Fagron	Isosorbide dinitrate ...	17C29-B01-350469	11
Fagron	Isosorbide dinitrate ...	17C29-B01-353161	1
Fagron	Isosorbide dinitrate ...	17C-B01-347589	1
Fagron	Isosorbide dinitrate ...	18J18-B01359579	1
Fagron	Isosorbide dinitrate ...	18J18-B01-358241	7
Fagron	Isosorbide dinitrate ...	17C29B01339447	1
Fagron	Isosorbide dinitrate ...	18J18B01358242	1
Fagron	Isosorbide dinitrate ...	18J18-B01	1
Fagron	Isosorbide dinitrate ...	17C29-B01-347589	9
Fagron	Isosorbide dinitrate ...	17C29B01347589	1
Fagron	Isosorbide dinitrate ...	17C29-B01-346333	1
Fagron	Isosorbide dinitrate ...	17C29-B01-347248	1
Fagron	Isosorbide dinitrate ...	17C29-b01-339448	1
Fagron	Isosorbide dinitrate ...	17C29-B01-352558	1
Fagron	Isosorbide dinitrate ...	17C29B01348859	1
Fagron	Isosorbide dinitrate ...	13A07-0003	1
Fagron	Isosorbide dinitrate ...	17C29B01335070	1
Fagron	Isosorbide dinitrate ...	101652	3
Fagron	Isosorbide dinitrate ...	12f08-n04	1
Fagron	Isosorbide dinitrate ...	14D09-B03_294282	1
Fagron	Isosorbide dinitrate ...	14D09-B03-29481	1
Fagron	Isosorbide dinitrate ...	14D09-B03-294282Eber	1
Fagron	Isosorbide dinitrate ...	14D09-B03-294282	3
Fagron	Isosorbide dinitrate ...	14D09-B03-300234	14
Fagron	Isosorbide dinitrate ...	17c29-b01-339448	1
Fagron	Isosorbide dinitrate ...	14-D09-B03-300234	1
Fagron	Isosorbide dinitrate ...	4314U-03516	1
Fagron	Isosorbide dinitrate ...	11F28-N05	2
Fagron	Isosorbide dinitrate ...	14D09-B03-308419	16
Fagron	Isosorbide dinitrate ...	8564340	1
Fagron	Isosorbide dinitrate ...	14D09B03300235	1
Fagron	Isosorbide dinitrate ...	14D09-B3-308419	1
Fagron	Isosorbide dinitrate ...	14D09-B03-313566	4
Fagron	Isosorbide dinitrate ...	1216M-03516	1
Fagron	Isosorbide dinitrate ...	14d09-b03-313566	1
Fagron	Isosorbide dinitrate ...	7013190003	1
Fagron	Isosorbide dinitrate ...	15/10-B01-315217	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Isosorbide dinitrate ...	14D09-B03_308419	1
Fagron	Isosorbide dinitrate ...	15 10-B01-315217	1
Fagron	Isosorbide dinitrate ...	15I10-B01-318522	1
Fagron	Isosorbide dinitrate ...	15/10B01-315217	1
Fagron	Isosorbide dinitrate ...	1510-B01-315217	2
Fagron	Isosorbide dinitrate ...	12k14-n02	1
Fagron	Isosorbide dinitrate ...	10C19-N01	1
Fagron	Isosorbide dinitrate ...	3211U-03516	1
Fagron	Isosorbide dinitrate ...	14D09-B03	8
Fagron	Isosorbide dinitrate ...	11J17-N05	3
Fagron	Isosorbide dinitrate ...	11D12-N01	3
Fagron	Isosorbide dinitrate ...	12C08-N05	4
Fagron	Isosorbide dinitrate ...	1309210009	1
Fagron	Isosorbide dinitrate ...	11I17-N05	1
Fagron	Isosorbide dinitrate ...	11/17-N05	2
Fagron	Isosorbide dinitrate ...	9101205	1
Fagron	Isosorbide dinitrate ...	10029-N10	1
Fagron	Isosorbide dinitrate ...	P11008	1
Fagron	Isosorbide dinitrate ...	10D29-N10	1
Fagron	Isosorbide dinitrate ...	12F08-N04	9
Fagron	Isosorbide dinitrate ...	12K14-N02	6
Fagron	Isosorbide dinitrate ...	15/10-B01-320081	11
Fagron	Isosorbide dinitrate ...	13A07-N13	8
Fagron	Isosorbide dinitrate ...	11117-N05	1
Fagron	Isosorbide dinitrate ...	13d15-n02	1
Fagron	Isosorbide dinitrate ...	3413U03516	1
Fagron	Isosorbide dinitrate ...	11/17-No5	1
Fagron	Isosorbide dinitrate ...	9E05-N05	1
Fagron	Isosorbide dinitrate ...	391311-03516	1
Fagron	Isosorbide dinitrate ...	12K14N02	1
Fagron	Isosorbide dinitrate ...	13F27-N04	7
Fagron	Isosorbide dinitrate ...	13I26-N18	10
Fagron	Isosorbide dinitrate ...	2114M3516	1
Fagron	Isosorbide dinitrate ...	14D09-B03-294280	9
Fagron	Isosorbide dinitrate ...	X12010	1
Fagron	Isosorbide dinitrate ...	2413Q03516	1
Fagron	Isosorbide dinitrate ...	15I10-B01-315217	2
Fagron	Isosorbide dinitrate ...	testxx	1
Fagron	Isosorbide dinitrate ...	14D09-B03-305633	1
Fagron	Isosorbide dinitrate ...	15-10-b01-325475	1
Fagron	Isosorbide dinitrate ...	25021711	1
Fagron	Isosorbide dinitrate ...	1510-B01320081	1
Fagron	Isosorbide dinitrate ...	1510-b01-328132	1
Fagron	Isosorbide dinitrate ...	15I10-B01-328132	2
Fagron	Isosorbide dinitrate ...	15/10-B01-328577	6
Fagron	Isosorbide dinitrate ...	15/10-b01-331287	3
Fagron	Isosorbide dinitrate ...	15/10-B01-328132	2
Fagron	Isosorbide dinitrate ...	15I10-B01-328577	3
Fagron	Isosorbide dinitrate ...	15/10-b01-328132	3
Fagron	Isosorbide dinitrate ...	1510B01328577	1
Fagron	Isosorbide dinitrate ...	15/10-B01-331287	1
Fagron	Isosorbide dinitrate ...	15-10-b01-331287	5
Fagron	Isosorbide dinitrate ...	15/10-B01-331286	2
Fagron	Isosorbide dinitrate ...	17c29-b01-335071	6
Fagron	Isosorbide dinitrate ...	15110-B01-331286	1
Fagron	Isosorbide dinitrate ...	1510B01325474	1
Fagron	Isosorbide dinitrate ...	17C29-B01	4
Fagron	Isosorbide dinitrate ...	193469	5
Fagron	Isosorbide dinitrate ...	17C29-B01335070	1
Fagron	Isosorbide dinitrate ...	17C29-B01-335071	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Isosorbide dinitrate ...	517I-03516	1
Fagron	Isosorbide dinitrate ...	15110-B01-325474	1
Fagron	Isosorbide dinitrate ...	151-b01-328132	2
Fagron	Isosorbide dinitrate ...	1510-B01-320081	1
Fagron	Isosorbide dinitrate ...	3115M-03516	1
Fagron	Isosorbide dinitrate ...	15-10-b01-318522	1
Fagron	Isosorbide dinitrate ...	160629I	1
Fagron	Isosorbide dinitrate ...	1510B01320081	1
Fagron	Isosorbide dinitrate ...	15I10-B01	8
Fagron	Isosorbide dinitrate ...	149168	10
Fagron	Isosorbide dinitrate ...	15I10-B01-320081	4
Fagron	Isosorbide dinitrate ...	13F27-N040234/1300	1
Fagron	Isosorbide dinitrate ...	15/10-B01-325474	1
Fagron	Isosorbide dinitrate ...	1510-b01-325475	3
Fagron	Isosorbide dinitrate ...	15I10-B01-310081	1
Fagron	Isosorbide dinitrate ...	15I10-B01-320782	1
Fagron	Isosorbide dinitrate ...	15I10-B01-320081	1
Fagron	Isosorbide dinitrate ...	15I10-B01-325474	5
Fagron	Isosorbide dinitrate ...	15110-B01-320081	1
Fagron	Ketamine hydrochloride	173017P	1
Fagron	Ketamine hydrochloride	14I03-B02-303165	1
Fagron	Ketamine hydrochloride	182513P	1
Fagron	Ketamine hydrochloride	14/03-B02-303165	1
Fagron	Ketamine hydrochloride	16E18-B03-328264	1
Fagron	Ketamine hydrochloride	14I03-B02-309852	1
Fagron	Ketamine hydrochloride	16E18-B03-332957	1
Fagron	Ketamine hydrochloride	16E18-B03	1
Fagron	Ketamine hydrochloride	16E18-B03-339897	2
Fagron	Ketamine hydrochloride	16E18-B03-334370	1
Fagron	Ketamine hydrochloride	15E07-B03-322164	5
Fagron	Ketamine hydrochloride	17K14-B01-358396	1
Fagron	Ketamine hydrochloride	17K14-B01-356831	1
Fagron	Ketamine hydrochloride	17K14-B01-347995	1
Fagron	Ketamine hydrochloride	17K14-B01-3533344	1
Fagron	Ketoconazole	17H14-B01-339422	1
Fagron	Ketoconazole	16A26-B01-335467	1
Fagron	Ketoconazole	12I17-N02	1
Fagron	Ketoconazole	16A26-B01-327490	1
Fagron	Ketoconazole	15G23-B03-316596	1
Fagron	Ketoconazole	15G23-B03-312063	1
Fagron	Ketoconazole	13H09-N07	1
Fagron	Ketoconazole	5015E-02866	1
Fagron	Ketoconazole	17H14-B01-344222	1
Fagron	Ketoconazole	14E09-B01-297637	1
Fagron	Ketoconazole	14E09-B01-297637	1
Fagron	Ketoconazole	14E09-B01-293736	1
Fagron	Ketoconazole	14E09-B01-293736	2
Fagron	Ketoconazole	13179613	1
Fagron	Ketoconazole	14K04-B06-309671	2
Fagron	Ketoconazole	13J23-N04	1
Fagron	Ketoconazole	17H14-B01-353376	1
Fagron	Lactose monohydrate	12C26-N05	1
Fagron	Lactose monohydrate	15A19-B06-303437	3
Fagron	Lactose monohydrate	14H27-B02	1
Fagron	Lactose monohydrate	14140312	14
Fagron	Lactose monohydrate	14A30-B07-296182	2
Fagron	Lactose monohydrate	14A30-B07	2
Fagron	Lactose monohydrate	13J23-B03-291114	1
Fagron	Lactose monohydrate	13I23-B03-287414	1
Fagron	Lactose monohydrate	13E23-B01-283334	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lactose monohydrate	13E23-B01	2
Fagron	Lactose monohydrate	1303C001-280388	1
Fagron	Lactose monohydrate	13E23-B01-283332	2
Fagron	Lactose monohydrate	282567/0	1
Fagron	Lactose monohydrate	13I23-B03	4
Fagron	Lactose monohydrate	11C29-N06	1
Fagron	Lactose monohydrate	12/23-B02-287554	1
Fagron	Lactose monohydrate	1301C006-278128	6
Fagron	Lactose monohydrate	1303C002-280398	1
Fagron	Lactose monohydrate	13/23-B02-287554	1
Fagron	Lactose monohydrate	1303C003-280390	1
Fagron	Lactose monohydrate	1303C001	1
Fagron	Lactose monohydrate	1301C005-278129	2
Fagron	Lactose monohydrate	1303C003	1
Fagron	Lactose monohydrate	11L22-N03	1
Fagron	Lactose monohydrate	13I23-B03-291116	1
Fagron	Lactose monohydrate	12G27-N02	5
Fagron	Lactose monohydrate	14H27-B01-305666	1
Fagron	Lactose monohydrate	12C09-N04	1
Fagron	Lactose monohydrate	1301C007-278131	1
Fagron	Lactose monohydrate	10726158	1
Fagron	Lactose monohydrate	10751120	1
Fagron	Lactose monohydrate	13I23-B03-291115	1
Fagron	Lactose monohydrate	14333525	1
Fagron	Lactose monohydrate	15A19-B06-303439	3
Fagron	Lactose monohydrate	14H27-B02-302033	1
Fagron	Lactose monohydrate	15H14-B03-314451	3
Fagron	Lactose monohydrate	1623123	1
Fagron	Lactose monohydrate	15A19-B06	6
Fagron	Ladys bedstraw herb,...	12369004	2
Fagron	Lavender flower	1489A120919	6
Fagron	Levocarnitine	9041504	1
Fagron	Levocarnitine	13111701	1
Fagron	Levocarnitine	13G19-N12	2
Fagron	Levothyroxin sodium	16G25-B99-324827	17
Fagron	L-glutathione	14H21-B02	1
Fagron	Lidocaine	LDB/116049	1
Fagron	Lidocaine	17G23-B04-350541	1
Fagron	Lidocaine	17K23B04353428	1
Fagron	Lidocaine	13J30-b02-287639	1
Fagron	Lidocaine	13J30-B02-287640	2
Fagron	Lidocaine	17K23-B04-350543	4
Fagron	Lidocaine	13J30-B02-287641	1
Fagron	Lidocaine	13J30-B02-292219	1
Fagron	Lidocaine	13J30-B02	4
Fagron	Lidocaine	5122014A	1
Fagron	Lidocaine	13j30B02292219	1
Fagron	Lidocaine	1014I-07006	1
Fagron	Lidocaine	17K23-B04-353428	1
Fagron	Lidocaine	14E13-B05-295549	1
Fagron	Lidocaine	14E13-B05	2
Fagron	Lidocaine	13E28-B40-289433	1
Fagron	Lidocaine	14E13-B05-295550	3
Fagron	Lidocaine	13H05-N05	1
Fagron	Lidocaine	17K23-B04-350541	3
Fagron	Lidocaine	17K23-B04-350189	3
Fagron	Lidocaine	13J30-B02-292219	1
Fagron	Lidocaine	LDB/117093	1
Fagron	Lidocaine	18J23-B02356573	1
Fagron	Lidocaine	13J30-B02-287641	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lidocaine	16F21-B06-323669	5
Fagron	Lidocaine	12E30-N10	10
Fagron	Lidocaine	11D28-N12	1
Fagron	Lidocaine	12C08-N11	1
Fagron	Lidocaine	11K17-N03	2
Fagron	Lidocaine	6J23-No2	1
Fagron	Lidocaine	12H16-N05	8
Fagron	Lidocaine	12I21-N06	4
Fagron	Lidocaine	12h16-N05	1
Fagron	Lidocaine		3
Fagron	Lidocaine	12H16-n05	1
Fagron	Lidocaine	12K13-N01	1
Fagron	Lidocaine	12/21-N06	2
Fagron	Lidocaine	13B27-N04	9
Fagron	Lidocaine	11C09-N07	1
Fagron	Lidocaine	13H05-N05	12
Fagron	Lidocaine	12h16_n05	1
Fagron	Lidocaine	214U-03584	1
Fagron	Lidocaine	13H05-N050171/245	1
Fagron	Lidocaine	1113M-03584	1
Fagron	Lidocaine	12I21N06	1
Fagron	Lidocaine	17K23-B04-347473	3
Fagron	Lidocaine	17K23B04347473	1
Fagron	Lidocaine	17G03-B08-341306	1
Fagron	Lidocaine	15H27-B07	1
Fagron	Lidocaine	15D23-B08	2
Fagron	Lidocaine	13B27-N04	1
Fagron	Lidocaine	16F21-B06-333787	1
Fagron	Lidocaine	15D23-B08_312058	1
Fagron	Lidocaine	15D23B08311708	1
Fagron	Lidocaine	15H27-B07-316784	1
Fagron	Lidocaine	15d23-b08-312058	2
Fagron	Lidocaine	16F21	2
Fagron	Lidocaine	16F21-B06-332253	1
Fagron	Lidocaine	15A15-B01	1
Fagron	Lidocaine	16F21B06332251	1
Fagron	Lidocaine	15H27-B07-315483	2
Fagron	Lidocaine	15h27-b05-318310	1
Fagron	Lidocaine	16A28-B06-320059	3
Fagron	Lidocaine	5742	1
Fagron	Lidocaine	16F21-B06-327232	2
Fagron	Lidocaine	16F21-B06	2
Fagron	Lidocaine	15D23-B08-312058	1
Fagron	Lidocaine	16F21B06325934	1
Fagron	Lidocaine	15D23-B08_309882	1
Fagron	Lidocaine	16F21-B06-325934	1
Fagron	Lidocaine	16F21-B06-328421	1
Fagron	Lidocaine	16A28B06320059	1
Fagron	Lidocaine	515M-07006	1
Fagron	Lidocaine	16F21-B06-333788	1
Fagron	Lidocaine	17K23-B04	5
Fagron	Lidocaine	17K23-B04-347474	2
Fagron	Lidocaine	17K23-B04-345215	3
Fagron	Lidocaine	13J30-B02-292220	1
Fagron	Lidocaine	14E13-B05-295549	2
Fagron	Lidocaine	1116E-07006	1
Fagron	Lidocaine	17K23-B04-343250	1
Fagron	Lidocaine	17G03-B08-343251	1
Fagron	Lidocaine	17G03-B08-341308	2
Fagron	Lidocaine	71-18	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lidocaine	15D23-B08-31058	1
Fagron	Lidocaine	17G03-B08-338286	1
Fagron	Lidocaine	17G03-B08	4
Fagron	Lidocaine	15A15-B01-304283	2
Fagron	Lidocaine	13J30-B02-292221	1
Fagron	Lidocaine	16f21-b06-335659	1
Fagron	Lidocaine	15A15B01304284	1
Fagron	Lidocaine	11C09-N07	1
Fagron	Lidocaine	15d23-b08-309882	1
Fagron	Lidocaine	15D23-B08-312058	9
Fagron	Lidocaine	15D23-B08-309882	1
Fagron	Lidocaine	16F21-B06-332253	3
Fagron	Lidocaine	15A15-B01-304282	4
Fagron	Lidocaine	17K23-B04-347475	7
Fagron	Lidocaine	16F21-B06-323670	15
Fagron	Lidocaine hydrochlori...	16C09-B03-319230	1
Fagron	Lidocaine hydrochlori...	16J24-B15-329421	1
Fagron	Lidocaine hydrochlori...	9B3M1L2Fo118	1
Fagron	Lidocaine hydrochlori...	17A16-B07-333010	4
Fagron	Lidocaine hydrochlori...	14e13-b06-294973	1
Fagron	Lidocaine hydrochlori...	11H18-N01	1
Fagron	Lidocaine hydrochlori...	12E14-N04	2
Fagron	Lidocaine hydrochlori...	LDH/116037	2
Fagron	Lidocaine hydrochlori...	18B23-B03-352117	3
Fagron	Lidocaine hydrochlori...	17A16-B07-331959	2
Fagron	Lidocaine hydrochlori...	218I-01775	1
Fagron	Lidocaine hydrochlori...	17G10-B08-342287	1
Fagron	Lidocaine hydrochlori...	13F07-N06	7
Fagron	Lidocaine hydrochlori...	18B23-B03-348590	2
Fagron	Lidocaine hydrochlori...	12H09-N08	2
Fagron	Lidocaine hydrochlori...	17G10-B08-346659	1
Fagron	Lidocaine hydrochlori...	168020-0001	1
Fagron	Lidocaine hydrochlori...	17G10-B08-346659	2
Fagron	Lidocaine hydrochlori...	13218-B11	1
Fagron	Lidocaine hydrochlori...	LDH/117022	1
Fagron	Lidocaine hydrochlori...	13F07N06	1
Fagron	Lidocaine hydrochlori...	14j20-b05-303133	1
Fagron	Lidocaine hydrochlori...	14J20-B05-303131	2
Fagron	Lidocaine hydrochlori...	14J07-B04-299433	3
Fagron	Lidocaine hydrochlori...	4612Q-01775	1
Fagron	Lidocaine hydrochlori...	15/17-B01316805	1
Fagron	Lidocaine hydrochlori...	12c05-N02	1
Fagron	Lidocaine hydrochlori...	14B12-B03-292214	11
Fagron	Lidocaine hydrochlori...	11L09-N09	2
Fagron	Lidocaine hydrochlori...	9E08-N03	1
Fagron	Lidocaine hydrochlori...	8F26-N07	1
Fagron	Lidocaine hydrochlori...	14E13-B06-294975	5
Fagron	Lidocaine hydrochlori...	13L18-B11-286901	1
Fagron	Lidocaine hydrochlori...	18B23-B03	2
Fagron	Lidocaine hydrochlori...	15I17-B01-315828	2
Fagron	Lidocaine hydrochlori...	15/17-B01-316805	7
Fagron	Lidocaine hydrochlori...	15I17-B01-316805	2
Fagron	Lidocaine hydrochlori...	16C09-B03	1
Fagron	Lidocaine hydrochlori...	14e13-b06-294975	1
Fagron	L-isoleucine	12080167	1
Fagron	L-isoleucine	13E02-B01-286930	1
Fagron	L-isoleucine	14D08-B02-296089	2
Fagron	L-isoleucine	17H14-B04-339008	3
Fagron	L-isoleucine	13A22-N03	1
Fagron	L-isoleucine	16H24-B03-331739	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Lithium orotate monoh...	13J04-N01	1
Fagron	Lungwort herb, cut	Nachprüfung	1
Fagron	Lysine hydrochloride	14H20-B03-301894	1
Fagron	Macrogols 4000	25082014A	1
Fagron	Macrogols 4000	15092014A	1
Fagron	Macrogols 4000	13L30-B31	1
Fagron	Macrogols 4000	20092014A	1
Fagron	Macrogols 4000	14A03B01	1
Fagron	Macrogols 4000	290714c	1
Fagron	Macrogols 4000	18082014C	1
Fagron	Macrogols 4000	13I30-N03	1
Fagron	Macrogols 4000	18G30-B01	3
Fagron	Macrogols 4000	14D14-B06	1
Fagron	Macrogols 4000	914I-02103	1
Fagron	Macrogols 4000	28052014B	1
Fagron	Macrogols 4000	22122014A	1
Fagron	Macrogols 4000	14A03-B01-289605	1
Fagron	Macrogols 4000	14J03-B07300840	1
Fagron	Macrogols 4000	14D14-B06-293453	2
Fagron	Macrogols 4000	14L12-B06	1
Fagron	Macrogols 4000	15H10-B04-314037	2
Fagron	Macrogols 4000	15G15-B07	1
Fagron	Macrogols 4000	15H10-B04-31407	1
Fagron	Macrogols 4000	15A08-B03-302761	1
Fagron	Macrogols 4000	15/29-B01-318038	2
Fagron	Macrogols 4000	15A08B03302762	1
Fagron	Macrogols 4000	2215I-02103	1
Fagron	Macrogols 4000	16A20-B15-320652	6
Fagron	Macrogols 4000	15I29-B01-318038	1
Fagron	Macrogols 4000	15376510	1
Fagron	Macrogols 4000	16a20-b15-320652	8
Fagron	Macrogols 4000	16A20-B15-320060	1
Fagron	Macrogols 4000	16152002	2
Fagron	Macrogols 4000	15G15-B07-313592	2
Fagron	Macrogols 4000	15A08-B03-302762	3
Fagron	Macrogols 4000	15B13-B04-305461	3
Fagron	Macrogols 4000	15b03-b08-304701	1
Fagron	Macrogols 4000	2022015A	1
Fagron	Macrogols 4000	10022015A	1
Fagron	Macrogols 4000	10022015B	1
Fagron	Macrogols 4000	14G21-B04-296899	1
Fagron	Macrogols 4000	14K03-B07-300522	3
Fagron	Macrogols 4000	14G15-B07-295632	1
Fagron	Macrogols 4000	505201C	1
Fagron	Macrogols 4000	14J03-B07-300840	1
Fagron	Macrogols 4000	21052015A	1
Fagron	Macrogols 4000	3114I-02103	1
Fagron	Macrogols 4000	8062015B	1
Fagron	Macrogols 4000	8062015C	1
Fagron	Macrogols 4000	22062015A	1
Fagron	Macrogols 4000	14K03-_b07-300522	1
Fagron	Macrogols 4000	13G19-N14	1
Fagron	Macrogols 4000	8012015A	1
Fagron	Macrogols 4000	26052014A	1
Fagron	Macrogols 4000	30092014B	1
Fagron	Macrogols 4000	16K25-B02-330527	5
Fagron	Macrogols 4000	13A10-N01	1
Fagron	Macrogols 4000	26082013A	1
Fagron	Macrogols 4000	13F12-N09	1
Fagron	Macrogols 4000	280913C	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Macrogols 4000	16102013A	1
Fagron	Macrogols 4000	13B20N01	1
Fagron	Macrogols 4000	2122013a	1
Fagron	Macrogols 4000	12F29-N01	2
Fagron	Macrogols 4000	2122013b	1
Fagron	Macrogols 4000	20122013C	1
Fagron	Macrogols 4000	17012014A	1
Fagron	Macrogols 4000	13F12N09	1
Fagron	Macrogols 4000	1112M-02103	1
Fagron	Macrogols 4000	3032014B	1
Fagron	Macrogols 4000	13I17-N03	1
Fagron	Macrogols 4000	14A03-B01-289602	3
Fagron	Macrogols 4000	13G19-N08	2
Fagron	Macrogols 4000	12L10-N06	2
Fagron	Macrogols 4000	13E22-N03	6
Fagron	Macrogols 4000	1803029-01	2
Fagron	Macrogols 4000	12I20-n03	1
Fagron	Macrogols 4000	12F13-N01	2
Fagron	Macrogols 4000	13130-N03	1
Fagron	Macrogols 4000	12F29-N01	1
Fagron	Macrogols 4000	160514a	1
Fagron	Macrogols 4000	13B20-N01	2
Fagron	Macrogols 6000 powder	12J22-N02	1
Fagron	Macrogols 6000 powder	11G14-N03	1
Fagron	Macrogols 6000 powder	14J17-B01	1
Fagron	Magnesium aspartate d...	17C21-B05-333111	1
Fagron	Magnesium carbonate, ...	12F20-N04	1
Fagron	Magnesium chloride he...	1301C10	3
Fagron	Magnesium chloride he...	12D11-N01	1
Fagron	Magnesium chloride he...	12h21-N02	1
Fagron	Magnesium chloride he...	1304C002	1
Fagron	Magnesium chloride he...	1304C001-281211	2
Fagron	Magnesium chloride he...	1301c009-277802	1
Fagron	Magnesium chloride he...	M1430	1
Fagron	Magnesium chloride he...	11K17-N10	1
Fagron	Magnesium chloride he...	1301C10-277801	1
Fagron	Magnesium chloride he...	26061305	1
Fagron	Magnesium chloride he...	12H21-N02	9
Fagron	Magnesium chloride he...	1304c004-281210	1
Fagron	Magnesium chloride he...	1304C003	1
Fagron	Magnesium chloride he...	1304C004-281210	1
Fagron	Magnesium chloride he...	1301C009-277802	4
Fagron	Magnesium chloride he...	12B16-N01	1
Fagron	Magnesium chloride he...	1301C011-277803	1
Fagron	Magnesium chloride he...	12H21-No2	1
Fagron	Magnesium chloride he...	M1625	1
Fagron	Magnesium chloride he...	13B14-B01-287292	1
Fagron	Magnesium chloride he...	12H21-N02	1
Fagron	Magnesium chloride he...	13411207	1
Fagron	Magnesium chloride he...	1304C003-281209	1
Fagron	Magnesium chloride he...	1301C012-277804	1
Fagron	Magnesium chloride he...	M1555	1
Fagron	Magnesium hydrogen ci...	16E02-B03-0326874	1
Fagron	Magnesium hydrogen ci...	14I12-B02	1
Fagron	Magnesium hydrogen ci...	13K11-B02-288161	1
Fagron	Magnesium hydrogen ci...	13K11-B02-...	1
Fagron	Magnesium hydrogen ci...	1056239	1
Fagron	Magnesium hydrogen ci...	13B22-N03	1
Fagron	Magnesium hydrogen ci...	13D30-N07	2
Fagron	Magnesium hydrogen ci...	16E02-B03-326874	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Magnesium hydrogen ci...	12G10-N01	2
Fagron	Magnesium hydrogen ci...	13K11-B02-288161	1
Fagron	Magnesium hydrogen ci...	17G24-B03-341019	1
Fagron	Magnesium hydrogen ci...	11H05-N07	1
Fagron	Magnesium hydrogen ci...	3047166	1
Fagron	Magnesium orotate di...	14E20-B05-296005	1
Fagron	Magnesium orotate di...	1000060461	1
Fagron	Magnesium orotate di...	13H06-N07	1
Fagron	Magnesium orotate di...	12G05-N03	1
Fagron	Magnesium orotate di...	16K22B01329783	1
Fagron	Magnesium orotate di...	160801701	1
Fagron	Magnesium orotate di...	14/13-B05-305616	1
Fagron	Magnesium stearate	12D26-B03-289435	1
Fagron	Magnesium sulfate hep...	12L10-N07	1
Fagron	Magnesium sulfate hep...	26251	1
Fagron	Magnesium sulfate hep...	31161	1
Fagron	Magnesium sulfate hep...	12G24-N02	1
Fagron	Magnesium sulfate hep...	12L10-N07	8
Fagron	Magnesium sulfate hep...	13H05-N02	1
Fagron	Mannitol	17I13-B11-341331	1
Fagron	Mannitol	17A27-B02-331292	4
Fagron	Mannitol	18J25-B08-356984	1
Fagron	Mannitol	17I13-B10	1
Fagron	Mannitol	18L03-B06-357782	1
Fagron	Mannitol	19A14B04360505	1
Fagron	Mannitol	17J30309-343413	1
Fagron	Mannitol	17/13-B04-340968	1
Fagron	Mannitol	17K29-B05	1
Fagron	Mannitol	720014-0001	1
Fagron	Mannitol	14G04-B07	2
Fagron	Mannitol	16E11-B07-322830	2
Fagron	Mannitol	16E11-B04	1
Fagron	Mannitol	17F27-B07	2
Fagron	Mannitol	18E25-B07	1
Fagron	Mannitol	16F27-B05-324516	1
Fagron	Mannitol	18G04-B13-354621	3
Fagron	Mannitol	17D18-B03-335437	2
Fagron	Mannitol	13K20-B01-287698	2
Fagron	Mannitol	17J30-B10-343414	1
Fagron	Mannitol	15K17-B01-315560	2
Fagron	Mannitol	15I24-B04	1
Fagron	Mannitol	17D18-B01	1
Fagron	Mannitol	17F27-B05-337913	2
Fagron	Mannitol	17D12-B01-334435	1
Fagron	Mannitol	18J04-B05-355467	1
Fagron	Mannitol	17F27B04337912	1
Fagron	Mannitol	14G04-B07-301692	10
Fagron	Mannitol	17F27-B07-338598	2
Fagron	Mannitol	16D20-B03-320486	1
Fagron	Mannitol	17I13-B04-340968	1
Fagron	Mannitol	17J23-B05-342007	2
Fagron	Mannitol	17/13-B11-341331	1
Fagron	Mannitol	17F27B04-337912	1
Fagron	Mannitol	17I13-B03-340967	2
Fagron	Mannitol	17J30-B09-343413	2
Fagron	Mannitol	17A12-B03	1
Fagron	Mannitol	14G04-B07-301693	1
Fagron	Mannitol	296196	1
Fagron	Mannitol	18L21-B04-358870	1
Fagron	Mannitol	14G004-B07-301692	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Mannitol	15B05-B04-304084	3
Fagron	Mannitol	18E25-B09-352755	2
Fagron	Mannitol	17A27-B02	1
Fagron	Mannitol	170612185	1
Fagron	Mannitol	17A12-B02-332068	1
Fagron	Mannitol	16/28-B01-328593	1
Fagron	Mannitol	17C28-B06	2
Fagron	Mannitol	18A30-B08-346849	1
Fagron	Mannitol	11219122	2
Fagron	Mannitol	15A26-B01-304717	5
Fagron	Mannitol	12K28-N01	1
Fagron	Mannitol	14G04-B07-296196	1
Fagron	Mannitol	1302C004	2
Fagron	Mannitol	16E11B07322830	1
Fagron	Mannitol	16H29-B02-326640	2
Fagron	Mannitol	17J04-B06-341798	1
Fagron	Mannitol	1302C005	1
Fagron	Mannitol	18B23-B10-349553	2
Fagron	Mannitol	13B01-B05-283812	1
Fagron	Mannitol	13K20-B01-294194	2
Fagron	Mannitol	13E24N01	1
Fagron	Mannitol	13K20-B01	1
Fagron	Mannitol	13K20-B01-294193	1
Fagron	Mannitol	16G27-B05-325659	1
Fagron	Mannitol	14G17-B06-301921	1
Fagron	Mannitol	17B08-B02	2
Fagron	Mannitol	16H25-B03-326639	1
Fagron	Mannitol	16L15-B04	1
Fagron	Mannitol	11L19-N15	1
Fagron	Mannitol	15K17-B03-314903	1
Fagron	Mannitol	720014-0002	1
Fagron	Mannitol	15K25-B06-318740	2
Fagron	Mannitol	17F27-B04-337912	1
Fagron	Mannitol	14G04-B07-301694	1
Fagron	Mannitol	16F27-B04-32483	1
Fagron	Mannitol	16H02-B01-325660	1
Fagron	Mannitol	16E11-B06	1
Fagron	Mannitol	16G27-B03-324517	1
Fagron	Mannitol	13K20-B01-287706	1
Fagron	Mannitol	17K29-B04-343940	1
Fagron	Mannitol	17I13B03340967	1
Fagron	Mannitol	16H29-B03-326641	1
Fagron	Mannitol	15K25-B04-317786	2
Fagron	Mannitol	16F01-B03-322189	1
Fagron	Mannitol	12E24-N03	1
Fagron	Mannitol	16B24-B09-320013	2
Fagron	Mannitol	18D09-B11-351436	2
Fagron	Mannitol	17A12-B02	1
Fagron	Matricaria flower		1
Fagron	Melatonin	13D12-N10	3
Fagron	Melatonin	16G06-B01-327131	3
Fagron	Melatonin	1305029-01	8
Fagron	Melatonin	13H29-B03-291180	1
Fagron	Melatonin	17I28-B04-344440	1
Fagron	Melatonin	17/28-B04-344440	2
Fagron	Melatonin	17C22-B01-337089	2
Fagron	Melatonin	11F09-N02	1
Fagron	Melatonin	13D12N10	1
Fagron	Melatonin	17C22-B01-337091	1
Fagron	Melatonin	13H21-N065x	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Melatonin	17C22-B01-337090	3
Fagron	Melatonin	16343702	1
Fagron	Melatonin	15C20-B01-306077	1
Fagron	Melatonin	17128-B04-344439	1
Fagron	Melatonin	13H29-B03	1
Fagron	Melatonin	15C20-B01-306078	3
Fagron	Melatonin	13e22-n12	1
Fagron	Melatonin	14G18-B01-297242	1
Fagron	Melatonin	17I28-B04-348747	1
Fagron	Melatonin	15G06-B03-315183	6
Fagron	Melatonin	13E22N12	1
Fagron	Melatonin	16G06-B01-329142	1
Fagron	Melatonin	15C20-B01-327131	1
Fagron	Melatonin	15C20-B01-308426	2
Fagron	Melatonin	14G18-B01-297244	1
Fagron	Melatonin	13E22-N12	1
Fagron	Melatonin	12K13-N02	1
Fagron	Melatonin	15G06-B03-318398	1
Fagron	Melatonin	17/28-B04-344438	2
Fagron	Melatonin	13H21-N06	3
Fagron	Melatonin	14G18-B01-297245	1
Fagron	Melatonin	13H29-B03-291180	1
Fagron	Melatonin	15C20-B01-308427	5
Fagron	Melatonin	15G06-B03-323733	1
Fagron	Melatonin	120009	1
Fagron	Menthol	1219-N10	1
Fagron	Menthol	13D16-N09	8
Fagron	Menthol	1509045-03	13
Fagron	Menthol	14100041-01	1
Fagron	Menthol	16L20T04044769	1
Fagron	Menthol	18L04-T06-65445	1
Fagron	Menthol	13016-N09	1
Fagron	Menthol	13D16-N090631/1300	1
Fagron	Menthol	13G09-N090204/520	1
Fagron	Menthol	13J21-N04	2
Fagron	Menthol	16L20-T04	1
Fagron	Menthol	976035	1
Fagron	Menthol	16L20-T04-053177	4
Fagron	Menthol	12C05-N03	1
Fagron	Menthol	778301	1
Fagron	Menthol	11C01-N01	1
Fagron	Menthol	14H29-M02	9
Fagron	Menthol	14H29-M04	1
Fagron	Menthol	12D23-N06	4
Fagron	Menthol	1807044-01	6
Fagron	Menthol	18G02-T30-062159	5
Fagron	Menthol	12F01-N14	2
Fagron	Menthol	12G20-N02	3
Fagron	Menthol	12I19-N10	2
Fagron	Menthol	12Cos-N03	1
Fagron	Menthol	11G19-N07	1
Fagron	Mepivacaine hydrochlo...	13G05N02	1
Fagron	Mepivacaine hydrochlo...	12F20-N07	1
Fagron	Mepivacaine hydrochlo...	12E31-B40-305570	1
Fagron	Mepivacaine hydrochlo...	15I16-B01	1
Fagron	Mesalazine	16I28-B02-327679	24
Fagron	Methacholine chloride	15F10-B02-315523	1
Fagron	Methacholine chloride	17H03-B02-346564	4
Fagron	Methacholine chloride	13G10-N05	1
Fagron	Methacholine chloride	17D13-B04-334378	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Methacholine chloride	17D13-B04	1
Fagron	Methacholine chloride	17D13-B04-334377	2
Fagron	Methacholine chloride	16L20-B01	1
Fagron	Methacholine chloride	16L20-B01-329970	1
Fagron	Methacholine chloride	15F10-B02-315831	1
Fagron	Methacholine chloride	13121712	1
Fagron	Methacholine chloride	312160	1
Fagron	Methacholine chloride	10031504	1
Fagron	Methacholine chloride	14K20-B03-300973	2
Fagron	Methacholine chloride	14G30-B03-296210	1
Fagron	Methacholine chloride	24061503	1
Fagron	Methacholine chloride	13/13-N03	2
Fagron	Methionine	15A15-B08-307753	1
Fagron	Methionine	13G09-N02	1
Fagron	Methoxsalen	12E02-N05	2
Fagron	Methoxsalen	15F23-B03	1
Fagron	Methoxsalen	15F23-B03-318166	5
Fagron	Methoxsalen	15F23-B03-309644	1
Fagron	Methoxsalen	14F06-B05-301978	1
Fagron	Methoxsalen	14F06-B05-296500	1
Fagron	Methoxsalen	14f06-b05-296500	1
Fagron	Methoxsalen	14F06-B05-296501	3
Fagron	Methoxsalen	15C10-B99-306266	1
Fagron	Methoxsalen	15F23-B03-316262	1
Fagron	Methoxsalen	MS01/13-14	1
Fagron	Methoxsalen	2111A-03262	2
Fagron	Methoxsalen	11F21-N02	1
Fagron	Methoxsalen	15J15-B03-321321	1
Fagron	Methoxsalen	12G02-N14	1
Fagron	Methoxsalen	11f21-n02	1
Fagron	Methoxsalen	21111814	1
Fagron	Methyl-4-hydroxybenzo...	20111125-1	1
Fagron	Methyl-4-hydroxybenzo...	12H02-N08	1
Fagron	Methyl-4-hydroxybenzo...	13B25-N03	2
Fagron	Methyl-4-hydroxybenzo...		1
Fagron	Methyl-4-hydroxybenzo...	12D18-N01	1
Fagron	Methyl-4-hydroxybenzo...	16256115	2
Fagron	Methyl-4-hydroxybenzo...	13E27-N09	1
Fagron	Methylcellulose	12H06-N06	2
Fagron	Methylcellulose	13E02-B03-290531	1
Fagron	Methylcellulose	15G08B14314845	1
Fagron	Methylcellulose	131301(09)	1
Fagron	Methylcellulose	13/26-N04	1
Fagron	Methylprednisolone	12D17-N04	1
Fagron	Methylprednisolone	18A03-B05	1
Fagron	Methylrosanilinium ch...	12I17-N05	1
Fagron	Methylrosanilinium ch...	12D19-N01	1
Fagron	Methylrosanilinium ch...	13H28-N03	18
Fagron	Methylrosanilinium ch...	13H28-No3	1
Fagron	Methylrosanilinium ch...	13H28-N03-0161-1030	1
Fagron	Methylrosanilinium ch...	13H28-N03	2
Fagron	Methylrosanilinium ch...	13G03-N011	1
Fagron	Methylrosanilinium ch...	12-17-N050109-140	1
Fagron	Methylrosanilinium ch...	13G03-N11	1
Fagron	Methylrosanilinium ch...	12H28_N06	1
Fagron	Methylrosanilinium ch...	13G03N11	1
Fagron	Methylrosanilinium ch...	1989	2
Fagron	Methylrosanilinium ch...	13H28-N03	1
Fagron	Metoprolol tartrate	17G03-B09-340197	1
Fagron	Metoprolol tartrate	1700305RII	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Metoprolol tartrate	16E31-B03-325535	1
Fagron	Metoprolol tartrate	18K26-B01-355911	1
Fagron	Metoprolol tartrate	17G03-B09-350669	1
Fagron	Metoprolol tartrate	17G03-B09-347089	1
Fagron	Metoprolol tartrate	16E31-B03-327107	23
Fagron	Metronidazole	14B27-B31-291234	1
Fagron	Metronidazole	14C14-B02-291241	1
Fagron	Metronidazole	13j21-N07	1
Fagron	Metronidazole	14C14-B02-291239	7
Fagron	Metronidazole	14C14-B02	4
Fagron	Metronidazole	48(17C10-B04-346277	1
Fagron	Metronidazole	14C14-B02-291237	9
Fagron	Metronidazole	14c14b02291276	1
Fagron	Metronidazole	14C14B02-291276	2
Fagron	Metronidazole	14c14-b02-291239	1
Fagron	Metronidazole	20140131f	1
Fagron	Metronidazole	14C14-Bo2-291276	1
Fagron	Metronidazole	14c14-b02-291275	1
Fagron	Metronidazole	14C14-B02-291276	5
Fagron	Metronidazole	17c10-b04-346276	1
Fagron	Metronidazole	17C10-B04-355865	1
Fagron	Metronidazole	172939-B01-295547	1
Fagron	Metronidazole	20140131F	3
Fagron	Metronidazole	14C14-B291237	1
Fagron	Metronidazole	14c14-b02-291276	2
Fagron	Metronidazole	17C10-B04-346279	1
Fagron	Metronidazole	14C14-B02-291238	2
Fagron	Metronidazole	14C14-802-291237	1
Fagron	Metronidazole	14C14-B02-296437	10
Fagron	Metronidazole	17C10-B04-346276	5
Fagron	Metronidazole	17C10-B04-346277	2
Fagron	Metronidazole	14C14-B02_291238	1
Fagron	Metronidazole	12f15-n04	1
Fagron	Metronidazole	14C14-B02-291275	4
Fagron	Metronidazole	12d30-n03	1
Fagron	Metronidazole	15a20-b07-308791	1
Fagron	Metronidazole	15A20-B07-312225	3
Fagron	Metronidazole	15A20-B07312057	1
Fagron	Metronidazole	16I19-B04-329159	2
Fagron	Metronidazole	1619B04329161	1
Fagron	Metronidazole	15A20-B07-312057	4
Fagron	Metronidazole	14C14-B01-295548	2
Fagron	Metronidazole	12A20-B07-303298	1
Fagron	Metronidazole	15A20/B07-303298	1
Fagron	Metronidazole	13I10-N03	3
Fagron	Metronidazole	14C14-B02-295603	1
Fagron	Metronidazole	15a20b07303298	1
Fagron	Metronidazole	15A20B07308525	1
Fagron	Metronidazole	141C14-B02-296437	1
Fagron	Metronidazole	15A20Bo7308196	1
Fagron	Metronidazole	14c14-B01-295545	1
Fagron	Metronidazole	14C14-B01	1
Fagron	Metronidazole	16/19-b04-329168	1
Fagron	Metronidazole	15a20-b07-308196	1
Fagron	Metronidazole	15A20-B07-303293	9
Fagron	Metronidazole	15A20-B07	1
Fagron	Metronidazole	12F15-N04	7
Fagron	Metronidazole	11F15-N07	2
Fagron	Metronidazole	14C14-B01-309058	2
Fagron	Metronidazole	606/10072017	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Metronidazole	15A20-B07-31225	1
Fagron	Metronidazole	15!20-B07-308791	1
Fagron	Metronidazole	16I19-B04-329167	3
Fagron	Metronidazole	160902	1
Fagron	Metronidazole	1619-B04-326155	1
Fagron	Metronidazole	16/19-B04-326154	1
Fagron	Metronidazole	17C10-B04-346275	2
Fagron	Metronidazole	16/19-B04	1
Fagron	Metronidazole	16J19-B04-326151	1
Fagron	Metronidazole	16/19bo4326155	1
Fagron	Metronidazole	16/16-B04-326151	1
Fagron	Metronidazole	16I19-B04-326151	4
Fagron	Metronidazole	16I19-B04-326150	1
Fagron	Metronidazole	16I19-B04	3
Fagron	Metronidazole	19I19-B04-326154	1
Fagron	Metronidazole	16/19-B04-329167	1
Fagron	Metronidazole	16I19-B04-326155	3
Fagron	Metronidazole	1619b04326155	1
Fagron	Metronidazole	16I19-B04-326150	3
Fagron	Metronidazole	16/19-B04-326150	7
Fagron	Metronidazole	16/19-B04-326151	2
Fagron	Metronidazole	16/19-B04-326152	1
Fagron	Metronidazole	17/19-304-329168	1
Fagron	Metronidazole	16/19-B04-329168	3
Fagron	Metronidazole	16D15-B50-324137	2
Fagron	Metronidazole	12a2o-n18	1
Fagron	Metronidazole	13/21-N7	2
Fagron	Metronidazole	16I19-B04-326152	2
Fagron	Metronidazole	17058003	9
Fagron	Metronidazole	15A20-B07-303294Mg	1
Fagron	Metronidazole	17C10	1
Fagron	Metronidazole	20130906F	4
Fagron	Metronidazole	13J21-N07	36
Fagron	Metronidazole	12D30-N02	1
Fagron	Metronidazole	17C10-B04	3
Fagron	Metronidazole	14C14-B02-300182	3
Fagron	Metronidazole	112356	1
Fagron	Metronidazole	1619-B04-326150	1
Fagron	Metronidazole	15A20-B07-308196	3
Fagron	Metronidazole	16/19-B04-329160	2
Fagron	Metronidazole	PM41	1
Fagron	Metronidazole	13/21-N07	1
Fagron	Metronidazole	1414-B02-291237	1
Fagron	Metronidazole	13D10-N02	4
Fagron	Metronidazole	12L17_N16	1
Fagron	Metronidazole	14-C14-B02-296437	1
Fagron	Metronidazole	14C14B02	3
Fagron	Metronidazole	15A20-B07-303295	1
Fagron	Metronidazole	13j21-n07	2
Fagron	Metronidazole	14C14-B02-2911239	1
Fagron	Metronidazole	13J21N07	2
Fagron	Metronidazole	14B27-B31-291235	1
Fagron	Metronidazole	16I19-B04-329162	1
Fagron	Metronidazole	605021726	2
Fagron	Metronidazole	15A20-B07-303294	2
Fagron	Metronidazole	12D30-N03	4
Fagron	Metronidazole	15A20-B07-308791	3
Fagron	Metronidazole	12L17N11	1
Fagron	Metronidazole	16I19-B04-329165	1
Fagron	Metronidazole	3212E01870	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Metronidazole	604081626	2
Fagron	Metronidazole	12L17-N16	4
Fagron	Metronidazole	11J07-N07	1
Fagron	Metronidazole	15A20-B07-308525	7
Fagron	Metronidazole	15A20-B07-308198	1
Fagron	Metronidazole	1220141155F	2
Fagron	Metronidazole	14C14-B01-306363	5
Fagron	Metronidazole	13B06-N03	11
Fagron	Metronidazole	11C21-N01	1
Fagron	Metronidazole	12L-N11	1
Fagron	Metronidazole	10F23-N05	1
Fagron	Metronidazole	16I19-B04-329161	1
Fagron	Metronidazole	13B06N03	1
Fagron	Metronidazole	9041303	1
Fagron	Metronidazole	12L17-n11	1
Fagron	Metronidazole	13D08N22	1
Fagron	Miconazole nitrate	18B06-B03-345422	1
Fagron	Miconazole nitrate	20081201	1
Fagron	Miconazole nitrate	18b06-b03-345419	1
Fagron	Miconazole nitrate	18B06-B03-347788	2
Fagron	Miconazole nitrate	12A16-N02	2
Fagron	Miconazole nitrate	15i08-b07-315481	1
Fagron	Miconazole nitrate	12H10-N02	4
Fagron	Miconazole nitrate	18C23-B02	3
Fagron	Miconazole nitrate	18C23-B02-350961	1
Fagron	Miconazole nitrate	15I08-B07-318686	3
Fagron	Miconazole nitrate	18C23-B02-351310	4
Fagron	Miconazole nitrate	15I08-B07-314838	1
Fagron	Miconazole nitrate	18C23-B02-248857	1
Fagron	Miconazole nitrate	18C23-B02-352488	1
Fagron	Miconazole nitrate	18C23-B02-353243	2
Fagron	Miconazole nitrate	18c23-b02	1
Fagron	Miconazole nitrate	18C23-B02-355295	1
Fagron	Miconazole nitrate	15/08-b07-318687	1
Fagron	Miconazole nitrate	15/08-b07-314838	2
Fagron	Miconazole nitrate	160711M	1
Fagron	Miconazole nitrate	12/21-M03	1
Fagron	Miconazole nitrate	27061202	1
Fagron	Miconazole nitrate	14d25-b06-299790	1
Fagron	Miconazole nitrate	17c01-b02-340778	1
Fagron	Miconazole nitrate	12K02-N06	7
Fagron	Miconazole nitrate	1316Q-03451	1
Fagron	Miconazole nitrate	15a13-b07-308432	1
Fagron	Miconazole nitrate	16E04-B03-321354	3
Fagron	Miconazole nitrate	15I08-B07-31485	1
Fagron	Miconazole nitrate	15/08-b07-315492	2
Fagron	Miconazole nitrate	15I08-B07-315492	1
Fagron	Miconazole nitrate	16E04-B03-327362	9
Fagron	Miconazole nitrate	16E04-B03-322131	6
Fagron	Miconazole nitrate	2900/03/16	1
Fagron	Miconazole nitrate	2900/04/14	2
Fagron	Miconazole nitrate	15/08-B07-316911	1
Fagron	Miconazole nitrate	14B18-B02-291247	1
Fagron	Miconazole nitrate	15I08-B07-215492	1
Fagron	Miconazole nitrate	2900/15/13	1
Fagron	Miconazole nitrate	14D25-B06-294959	8
Fagron	Miconazole nitrate	16e04-b03-322131	2
Fagron	Miconazole nitrate	14D25-B06	4
Fagron	Miconazole nitrate	14B18-B02-291250	1
Fagron	Miconazole nitrate	18C23-B02-347802	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Miconazole nitrate	1508-B07-318686	1
Fagron	Miconazole nitrate	14B18B02	1
Fagron	Miconazole nitrate	16E04-B03	3
Fagron	Miconazole nitrate	14b18-b02-291247	1
Fagron	Miconazole nitrate	16E04-B03-322132	5
Fagron	Miconazole nitrate	15I08-B07-312321	7
Fagron	Miconazole nitrate	15A13-B07-307140	6
Fagron	Miconazole nitrate	15I08-B07	3
Fagron	Miconazole nitrate	14H06-B01-299798	1
Fagron	Miconazole nitrate	14b18-B02-291231	1
Fagron	Miconazole nitrate	15A13-B07-306791nron	1
Fagron	Miconazole nitrate	16E04-B03-322124	13
Fagron	Miconazole nitrate	15a13-b07-307140	1
Fagron	Miconazole nitrate	15A13-807	1
Fagron	Miconazole nitrate	15/08-B07-312317	1
Fagron	Miconazole nitrate	15A13-B07-306791	8
Fagron	Miconazole nitrate	15I08-B07-312321	1
Fagron	Miconazole nitrate	14H06-B01-299799	8
Fagron	Miconazole nitrate	2900/16/14	2
Fagron	Miconazole nitrate	14D25-B06-299790	2
Fagron	Miconazole nitrate	14H06-B01-299813	3
Fagron	Miconazole nitrate	14D25-B06-294958	11
Fagron	Miconazole nitrate	14H16-B01-299798	1
Fagron	Miconazole nitrate	14149313	1
Fagron	Miconazole nitrate	16E04-B03-330082	6
Fagron	Miconazole nitrate	1508-B07-314838	1
Fagron	Miconazole nitrate	14B18-B02	4
Fagron	Miconazole nitrate	15I08-b07-312321	1
Fagron	Miconazole nitrate	15/08-B07-318686	2
Fagron	Miconazole nitrate	17C01-B02-340778	3
Fagron	Miconazole nitrate	4111301	1
Fagron	Miconazole nitrate	31101306	1
Fagron	Miconazole nitrate	13H28-N06	13
Fagron	Miconazole nitrate	13H28-N06	2
Fagron	Miconazole nitrate	17c01-b02-334877	1
Fagron	Miconazole nitrate	13A14-B09-287453	1
Fagron	Miconazole nitrate	13E13-N04	2
Fagron	Miconazole nitrate	17C01-B02-341395	4
Fagron	Miconazole nitrate	13D22-N14	7
Fagron	Miconazole nitrate	13E13-N08	2
Fagron	Miconazole nitrate	18B06-B03-345419	5
Fagron	Miconazole nitrate	17011307	1
Fagron	Miconazole nitrate	18b06-b03-345423	1
Fagron	Miconazole nitrate	18B06-B03-345423	2
Fagron	Miconazole nitrate	17C01-B02	2
Fagron	Miconazole nitrate	13A14B09283005	1
Fagron	Miconazole nitrate	15I08-B07-318485	1
Fagron	Miconazole nitrate	15/08-B07-315481	2
Fagron	Miconazole nitrate	1328-N06	1
Fagron	Miconazole nitrate	14B18-B02-291232	4
Fagron	Miconazole nitrate	17C01-B02-334877	4
Fagron	Miconazole nitrate	14B18-B02-291249	1
Fagron	Miconazole nitrate	2900/09/16	3
Fagron	Miconazole nitrate	17C01-B02-225671	1
Fagron	Miconazole nitrate	13A14-B09-287488	1
Fagron	Miconazole nitrate	14B18-B02-291248	6
Fagron	Miconazole nitrate	16E04-B03-321367	3
Fagron	Miconazole nitrate	15I08-B07-3121321	1
Fagron	Miconazole nitrate	15i08-b07-312317	1
Fagron	Miconazole nitrate	2900/08/12	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Miconazole nitrate	13H28_N06	1
Fagron	Miconazole nitrate	13A14-B09-287487	3
Fagron	Miconazole nitrate	15108-B07-318485	1
Fagron	Miconazole nitrate	17C01-B02-335671	4
Fagron	Miconazole nitrate	15108-B07-315492	3
Fagron	Miconazole nitrate	13A14-B09-287486	5
Fagron	Miconazole nitrate	15108-B07-318686	1
Fagron	Miconazole nitrate	17c01-b02-335671	1
Fagron	Miconazole nitrate	13A14-B09	2
Fagron	Miconazole nitrate	16E04-B03-332373	4
Fagron	Minoxidil	18C27-B01-354524	1
Fagron	Minoxidil	18C20-B01-352629	2
Fagron	Minoxidil	18C27-B01-354525	2
Fagron	Minoxidil	18C20-b01-351844	1
Fagron	Minoxidil	17H08-B09-346945	1
Fagron	Minoxidil	18C20-B01-352702	1
Fagron	Minoxidil	15436408	4
Fagron	Minoxidil	17H08-B09-346948	1
Fagron	Minoxidil	15K17-B06-324965	2
Fagron	Minoxidil	12C29-N050087/240	1
Fagron	Minoxidil	15K17-B06-321927	2
Fagron	Minoxidil	15K17-B05-318411	2
Fagron	Minoxidil	1520002631	1
Fagron	Minoxidil	15K197-B05-318411	1
Fagron	Minoxidil	15K17-B05-318691	2
Fagron	Minoxidil	15K17-B05-320188	3
Fagron	Minoxidil	15K17-B05-318240	2
Fagron	Minoxidil	15H07-B01-314726	1
Fagron	Minoxidil	15K17-B05-318241	1
Fagron	Minoxidil	1510001251	3
Fagron	Minoxidil	15E19-B09-308181	1
Fagron	Minoxidil	15F02-B09-312663	2
Fagron	Minoxidil	1520000768	1
Fagron	Minoxidil	12H07-N030265/432	1
Fagron	Minoxidil	15E19-B06-308183	1
Fagron	Minoxidil	15K17-B06-325490	1
Fagron	Minoxidil	15K17-B06-321928	2
Fagron	Minoxidil	11E17-N07	1
Fagron	Minoxidil	612212844	1
Fagron	Minoxidil	12C29-N05	6
Fagron	Minoxidil	120303503	2
Fagron	Minoxidil	11H07-N03	1
Fagron	Minoxidil	12H07N03	1
Fagron	Minoxidil	13K04-B05-290380	3
Fagron	Minoxidil	17B22-B01-334447	1
Fagron	Minoxidil	16J24-B04-328218	2
Fagron	Minoxidil	15K17-B06-325491	2
Fagron	Minoxidil	12H07-N03	1
Fagron	Minoxidil	16/24-B04-328218	1
Fagron	Minoxidil	15K17-B06	1
Fagron	Minoxidil	13H07-N03	1
Fagron	Minoxidil	321926	1
Fagron	Minoxidil	15E19-B06-308179	1
Fagron	Minoxidil	12H07-N03	1
Fagron	Minoxidil	15E19-B06-308180	5
Fagron	Minoxidil	27081301	1
Fagron	Minoxidil	11021402	1
Fagron	Minoxidil	13J29-B02	1
Fagron	Minoxidil	13J23N07	2
Fagron	Minoxidil	14L03-Bo2-301430	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Minoxidil	13123-N07	1
Fagron	Minoxidil	13K04-B05-290377	5
Fagron	Minoxidil	15E19-B06	1
Fagron	Minoxidil	13K04-B05-290379	3
Fagron	Minoxidil	1420001962	1
Fagron	Minoxidil	14L03-B02-301429	5
Fagron	Minoxidil	14L03-B02-301430	2
Fagron	Minoxidil	14L03-B02-301541	1
Fagron	Minoxidil	14B04-B07-297364	1
Fagron	Minoxidil	13K04-B05-0004	1
Fagron	Minoxidil	14B04-B07-290883	3
Fagron	Minoxidil	15E19-B06-308182	1
Fagron	Minoxidil	14L03-B02301434	1
Fagron	Minoxidil	13I23-N07	8
Fagron	Minoxidil	14L03-B02-301431	1
Fagron	Minoxidil	15E19-06-308183	1
Fagron	Minoxidil	15e19-b06-308180	1
Fagron	Minoxidil	15E19-B06-308181	2
Fagron	Minoxidil	15E19-B06-308183	6
Fagron	Minoxidil	14L03-B02	2
Fagron	Minoxidil	14L03-B02-30129	1
Fagron	Minoxidil	1520000768	1
Fagron	Minoxidil	4110I-02844	1
Fagron	Minoxidil	15F02-B08	1
Fagron	Minoxidil	13E13N10	1
Fagron	Minoxidil	18C20-B01-352116	1
Fagron	Minoxidil	12H07-N03	1
Fagron	Minoxidil	14L03-B02301429	1
Fagron	Minoxidil	26111302	1
Fagron	Minoxidil	12C29-N050145/240	1
Fagron	Minoxidil	14L03-B02-301429	1
Fagron	Minoxidil	15E19-B06-308182	5
Fagron	Minoxidil	170612184	1
Fagron	Minoxidil	7/13	1
Fagron	Minoxidil	12A25-N04	4
Fagron	Minoxidil	17B22-B02-336689	2
Fagron	Minoxidil	17B22-B02	1
Fagron	Minoxidil	16J24-B04-328226	3
Fagron	Minoxidil	17B22-B02-338767	1
Fagron	Minoxidil	17H08-B06	2
Fagron	Minoxidil	16j24-B04-329291	1
Fagron	Minoxidil	17B22-B02-338766	1
Fagron	Minoxidil	17H08-B06-341376	5
Fagron	Minoxidil	17H08-B06-341377	5
Fagron	Minoxidil	16J24-B04-329291	1
Fagron	Minoxidil	18C19-B01-348985	1
Fagron	Minoxidil	18C20-B01	1
Fagron	Minoxidil	18C19-B01-348984	3
Fagron	Minoxidil	17H08-B06-340514	1
Fagron	Minoxidil	18C20-B01-351844	1
Fagron	Minoxidil	17B22-B02-336690	1
Fagron	Minoxidil	17B22-B01-335427	4
Fagron	Minoxidil	17B22/B02/336689	2
Fagron	Minoxidil	16J24-B05	1
Fagron	Minoxidil	2/16M	2
Fagron	Minoxidil	17B22-B01	2
Fagron	Minoxidil	1/16M	1
Fagron	Minoxidil	11K11-N09	3
Fagron	Mometasone furoate, m...	14C07-B05-291266	2
Fagron	Mometasone furoate, m...	15/11-B02-316705	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Mometasone furoate, m...	16A20-B06-319748	7
Fagron	Mometasone furoate, m...	15G08-B11-322737	4
Fagron	Mometasone furoate, m...	16A20-B06-322739	4
Fagron	Mometasone furoate, m...	15B11-B03-308860	2
Fagron	Mometasone furoate, m...	MOM008/0915	1
Fagron	Mometasone furoate, m...	NS1513M	1
Fagron	Mometasone furoate, m...	18B19-B07-351351	2
Fagron	Mometasone furoate, m...	17D19B02342642	1
Fagron	Mometasone furoate, m...	15I11-B02-315705	1
Fagron	Mometasone furoate, m...	15G08B11310864	1
Fagron	Mometasone furoate, m...	140305MF	1
Fagron	Mometasone furoate, m...	14E20-B04-291265	3
Fagron	Mometasone furoate, m...	14E20-B04-294645	5
Fagron	Mometasone furoate, m...	5/0715	1
Fagron	Mometasone furoate, m...	10E31-N26	1
Fagron	Mometasone furoate, m...	16J24-B11-332110	1
Fagron	Mometasone furoate, m...	15G08-B11-0002	1
Fagron	Mometasone furoate, m...	17D19-B02339277	1
Fagron	Mometasone furoate, m...	18B19-B07	3
Fagron	Mometasone furoate, m...	2114M3640	1
Fagron	Mometasone furoate, m...	14H06-B07297240	1
Fagron	Mometasone furoate, m...	14H06-B07-298482	7
Fagron	Mometasone furoate, m...	14H06-B07-302375	2
Fagron	Mometasone furoate, m...	14e20-b04-291265	1
Fagron	Mometasone furoate, m...	335793	1
Fagron	Mometasone furoate, m...	13j23-no2	1
Fagron	Mometasone furoate, m...	16J24-B11-332111	3
Fagron	Mometasone furoate, m...	3313A-07030	1
Fagron	Mometasone furoate, m...	NU0698M	3
Fagron	Mometasone furoate, m...	15G08-B11	2
Fagron	Mometasone furoate, m...	15J11-B02-315705	1
Fagron	Mometasone furoate, m...	2215I-07030	1
Fagron	Mometasone furoate, m...	17D19-B02-342643	19
Fagron	Mometasone furoate, m...	13Do8-N06	1
Fagron	Mometasone furoate, m...	15G08-B11-317182	9
Fagron	Mometasone furoate, m...	15G08-B11-322738	5
Fagron	Mometasone furoate, m...	11H10-N04	1
Fagron	Mometasone furoate, m...	NU2504M	1
Fagron	Mometasone furoate, m...	14E20-B04-291267	3
Fagron	Mometasone furoate, m...	16A20-B06	3
Fagron	Mometasone furoate, m...	15B11-B03-304120	6
Fagron	Mometasone furoate, m...	17D19-B02-347348	4
Fagron	Mometasone furoate, m...	15/11B02315705	2
Fagron	Mometasone furoate, m...	15B11B03304122	2
Fagron	Mometasone furoate, m...	347349	2
Fagron	Mometasone furoate, m...	15B11-B03-304121	7
Fagron	Mometasone furoate, m...	15B11-B03-304122	5
Fagron	Mometasone furoate, m...	14H06-B07-302376	2
Fagron	Mometasone furoate, m...	14E20-B04-302374	2
Fagron	Mometasone furoate, m...	17D19-B02-347349	5
Fagron	Mometasone furoate, m...	17d19.b02.347349	1
Fagron	Mometasone furoate, m...	14H06-B07-298481	3
Fagron	Mometasone furoate, m...	2714U-07030	1
Fagron	Mometasone furoate, m...	16A20-B06-319749	6
Fagron	Mometasone furoate, m...	18B19-B07-341351	1
Fagron	Mometasone furoate, m...	16J24-B11-335793	11
Fagron	Mometasone furoate, m...	18B19-B07-352632	6
Fagron	Mometasone furoate, m...	14H06-B07-297240	1
Fagron	Mometasone furoate, m...	14C07-B05-291264	1
Fagron	Mometasone furoate, m...	15B11-B03-304124	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Mometasone furoate, m...	15B11-B03-304123	12
Fagron	Mometasone furoate, m...	15B11-B03	2
Fagron	Mometasone furoate, m...	15/11-B02	4
Fagron	Mometasone furoate, m...	15I11-B02-314479	4
Fagron	Mometasone furoate, m...	16J24B11334371	1
Fagron	Mometasone furoate, m...	15B11B03308860	1
Fagron	Mometasone furoate, m...	14e20b04291265	1
Fagron	Mometasone furoate, m...	15/11-B02-315705	3
Fagron	Mometasone furoate, m...	16A20-B06319749	1
Fagron	Mometasone furoate, m...	16J24B11332112	1
Fagron	Mometasone furoate, m...	15B11B03304121	1
Fagron	Mometasone furoate, m...	334371	1
Fagron	Mometasone furoate, m...	15/11-B02-314479	6
Fagron	Mometasone furoate, m...	15G08-B11-312233	4
Fagron	Mometasone furoate, m...	15G-B11-310864	1
Fagron	Mometasone furoate, m...	15B11-B03-308276	1
Fagron	Mometasone furoate, m...	15G08-B11-310864	10
Fagron	Mometasone furoate, m...	17d19-b02-347349	1
Fagron	Mometasone furoate, m...	17D19_B02_339278	1
Fagron	Mometasone furoate, m...	15b11-b03-308276	1
Fagron	Mometasone furoate, m...	4(18B19-B07-351001	1
Fagron	Mometasone furoate, m...	18B19-B07-351001	1
Fagron	Mometasone furoate, m...	16J24-11-335793	1
Fagron	Mometasone furoate, m...	701834-0004	1
Fagron	Mometasone furoate, m...	13j23-N02	1
Fagron	Mometasone furoate, m...	16C24-B03-326790	7
Fagron	Mometasone furoate, m...	12B22-N01	1
Fagron	Mometasone furoate, m...	12B22-N04	1
Fagron	Mometasone furoate, m...	10C29-N10	1
Fagron	Mometasone furoate, m...	18819-b07	1
Fagron	Mometasone furoate, m...	13F17-N01	7
Fagron	Mometasone furoate, m...	17D19-B02	1
Fagron	Mometasone furoate, m...	11H10-N01	1
Fagron	Mometasone furoate, m...	16C24B03	1
Fagron	Mometasone furoate, m...	11K28-N09	1
Fagron	Mometasone furoate, m...	13H13-N04	1
Fagron	Mometasone furoate, m...	10805-N12	1
Fagron	Mometasone furoate, m...	16C24-B03-328272	5
Fagron	Mometasone furoate, m...	16J24-B11-332112	12
Fagron	Mometasone furoate, m...	13d08-n06	1
Fagron	Mometasone furoate, m...	12I29-N03	1
Fagron	Mometasone furoate, m...	16C24-B03	3
Fagron	Mometasone furoate, m...	13D08-N06	17
Fagron	Mometasone furoate, m...	16C24-B03-327346	4
Fagron	Mometasone furoate, m...	11h10-n04	1
Fagron	Mometasone furoate, m...	18B19-B07-352744	1
Fagron	Mometasone furoate, m...	16J24-B11-329817	4
Fagron	Mometasone furoate, m...	16C24B03324851	1
Fagron	Mometasone furoate, m...	18819-B07-357640	1
Fagron	Mometasone furoate, m...	12H09-N02	7
Fagron	Mometasone furoate, m...	16C24	1
Fagron	Mometasone furoate, m...	16J24-B11-334371	6
Fagron	Mometasone furoate, m...	16C24-b03-327372	1
Fagron	Mometasone furoate, m...	16/24-B11-332122	1
Fagron	Mometasone furoate, m...	13F7-N01	1
Fagron	Mometasone furoate, m...	12E02-N11	2
Fagron	Mometasone furoate, m...	352632	1
Fagron	Mometasone furoate, m...	16j24-b11-329817	1
Fagron	Mometasone furoate, m...	16j24-b11	1
Fagron	Mometasone furoate, m...	12J29-N03	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Mometasone furoate, m...	16C24-B03-327372	4
Fagron	Mometasone furoate, m...	12J29-N01	1
Fagron	Mometasone furoate, m...	18B19-B07-357640	1
Fagron	Mometasone furoate, m...	16J24-B11	5
Fagron	Mometasone furoate, m...	12e02-n11	1
Fagron	Mometasone furoate, m...	13D08N06	1
Fagron	Mometasone furoate, m...	18B19-B07-354320	5
Fagron	Mometasone furoate, m...	13f17-n01	1
Fagron	Mometasone furoate, m...	16C24-B03-324852	6
Fagron	Mometasone furoate, m...	16C24B03324852	1
Fagron	Mometasone furoate, m...	18B19-B07-352745	2
Fagron	Mometasone furoate, m...	15G08	1
Fagron	Mometasone furoate, m...	18B19B07352745	1
Fagron	Mometasone furoate, m...	16724-B11-332112	1
Fagron	Mometasone furoate, m...	16A20-B06-319740	1
Fagron	Mometasone furoate, m...	16C24-803-324851	1
Fagron	Mometasone furoate, m...	1814E3640	1
Fagron	Mometasone furoate, m...	11C02-N03	2
Fagron	Mometasone furoate, m...	11F08-N03	1
Fagron	Mometasone furoate, m...	16C24-B03-324851	5
Fagron	Mometasone furoate, m...	15G08B11322738	1
Fagron	Mupirocin calcium	15H13-B01	1
Fagron	Mupirocin calcium	14B03-B30-297291	1
Fagron	Mupirocin calcium	17F21-B03-338597	1
Fagron	Mupirocin calcium	83200340116	1
Fagron	Mupirocin calcium	17F21-B03-339080	1
Fagron	Naloxone hydrochlorid...	14F20-B02-296077	1
Fagron	Naloxone hydrochlorid...	13L12-B03-289511	1
Fagron	Naloxone hydrochlorid...	15B06-B02-304421	1
Fagron	Naloxone hydrochlorid...	15H04-B09-314068	1
Fagron	Naloxone hydrochlorid...	15H04-B09	3
Fagron	Naloxone hydrochlorid...	12G19-N02	1
Fagron	Naloxone hydrochlorid...	12-00071	2
Fagron	Naloxone hydrochlorid...	12G19-N02	2
Fagron	Naloxone hydrochlorid...	15E28-B01-308324	3
Fagron	Naloxone hydrochlorid...	11/23-N01	1
Fagron	Naloxone hydrochlorid...	15K11-B05-319571	1
Fagron	Naphazoline hydrochlo...	11K25-N03	1
Fagron	Naphazoline hydrochlo...	15J05-B05-322527	1
Fagron	Naphazoline hydrochlo...	13K27-N04	11
Fagron	Naphazoline hydrochlo...	1802201903	1
Fagron	Naphazoline hydrochlo...	18D09-B08-350950	3
Fagron	Naphazoline hydrochlo...	17C02-B05-346657	2
Fagron	Naphazoline hydrochlo...	17C02-B05-343566	4
Fagron	Naphazoline hydrochlo...	17C02-B05-341579	1
Fagron	Naphazoline hydrochlo...	7001260012	1
Fagron	Naphazoline hydrochlo...	17C02-B05-335113	1
Fagron	Naphazoline hydrochlo...	13I03-B03-305291	1
Fagron	Naphazoline hydrochlo...	17c02-b05-335313	1
Fagron	Naphazoline hydrochlo...	13/03-B03-305291	1
Fagron	Naphazoline hydrochlo...	13I03-B03-295105	2
Fagron	Naphazoline hydrochlo...	1303-B03	1
Fagron	Naphazoline hydrochlo...	13I03-B03-295105	2
Fagron	Naphazoline hydrochlo...	13/03-B03305291	1
Fagron	Naphazoline hydrochlo...	13C11-N07	1
Fagron	Naphazoline hydrochlo...	13k27-n04	1
Fagron	Naphazoline hydrochlo...	60494000113K27-N04	1
Fagron	Naphazoline hydrochlo...	15J05-B05-326850	3
Fagron	Naphazoline hydrochlo...	13K27-N04	1
Fagron	Naphazoline hydrochlo...	17C02-B05-335313	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Naphazoline hydrochlo...	15J05-B05-330834	3
Fagron	Naphazoline hydrochlo...	17C02-B05-332933	4
Fagron	Naphazoline hydrochlo...	13103-803-305291	1
Fagron	Naphazoline hydrochlo...	18D09B08352462	1
Fagron	Naphazoline hydrochlo...	18D09-B08-354882	1
Fagron	Naphazoline hydrochlo...	13A07-N06	4
Fagron	Naphazoline hydrochlo...	4309I-01885	1
Fagron	Naphazoline hydrochlo...	11J13-N08	1
Fagron	Naphazoline hydrochlo...	11K25-N03	1
Fagron	Naphazoline hydrochlo...	13103-B03-295105	1
Fagron	Naphazoline hydrochlo...	11/05-N02	1
Fagron	Naphazoline hydrochlo...	12K29-N01o256/310	1
Fagron	Naphazoline hydrochlo...	13C11N07	1
Fagron	Naphazoline hydrochlo...	13/03-B03	1
Fagron	Naphazoline hydrochlo...	13C11-N07	4
Fagron	Naphazoline hydrochlo...	13I03-B03	1
Fagron	Naphazoline hydrochlo...	15J05B05326850	1
Fagron	Naphazoline hydrochlo...	11I05-N02	2
Fagron	Naphazoline hydrochlo...	13I03-B03-305291	1
Fagron	Naphazoline hydrochlo...	13/03-B03-295105	1
Fagron	Naphazoline hydrochlo...	15J05-B05	2
Fagron	Naphazoline hydrochlo...	13I03-B03	1
Fagron	Naphazoline hydrochlo...	17C02-B05	3
Fagron	Naphazoline hydrochlo...	4509E-01885	2
Fagron	Naphazoline nitrate	16A14-B03-327343	1
Fagron	Naphazoline nitrate	12J08-B40	1
Fagron	Naphazoline nitrate	12J08B40-294166	1
Fagron	Naphazoline nitrate	11I05-N03	1
Fagron	Naphazoline nitrate	12/21-N05	1
Fagron	Naphazoline nitrate	16A14-B03-322534	1
Fagron	Naphazoline nitrate	12J08-B40-294166	1
Fagron	Naphazoline nitrate	140214N	1
Fagron	Naphazoline nitrate	21021304	1
Fagron	Naphazoline nitrate	12j08-b40-303401	1
Fagron	Naproxen	16B04-B09-321569	24
Fagron	Naproxen	13E31-B40-287536	2
Fagron	Naproxen	17C10-B03-336472	1
Fagron	Neomycin sulfate	18G28-B05-354472	1
Fagron	Neomycin sulfate	12H23-N04	1
Fagron	Neomycin sulfate	M71862	2
Fagron	Neomycin sulfate	F08469	1
Fagron	Neomycin sulfate	13K15-B33-288073	1
Fagron	Neomycin sulfate	13D10-N05	1
Fagron	Neomycin sulfate	16A20-B01-327433	1
Fagron	Neomycin sulfate	13G09N07	2
Fagron	Neomycin sulfate	16A20-B10-327433	2
Fagron	Neomycin sulfate	13K15-B33-288069	2
Fagron	Neomycin sulfate	13G09-N070129/340	1
Fagron	Neomycin sulfate	15J09-B04-314017	1
Fagron	Neomycin sulfate	15J09-B04	1
Fagron	Neomycin sulfate	13D10-N05	6
Fagron	Neomycin sulfate	13G09-N06	1
Fagron	Neomycin sulfate	15J09-B04314015	1
Fagron	Neomycin sulfate	13G09-N07	3
Fagron	Neomycin sulfate	11C21-N10	4
Fagron	Neomycin sulfate	15J27-B07	1
Fagron	Neomycin sulfate	14A31-B02	1
Fagron	Neomycin sulfate	15J09-B04-314015	4
Fagron	Neomycin sulfate	15J09-B04-314017	5
Fagron	Neomycin sulfate	14C27-B05-302517	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Neomycin sulfate	14C27-B05-302516	2
Fagron	Neomycin sulfate	14C27-B05	1
Fagron	Neomycin sulfate	NMZ11006K	1
Fagron	Neomycin sulfate	15i09-B04-314015	1
Fagron	Neomycin sulfate	12H23-N04	5
Fagron	Neomycin sulfate	11/19-N14	1
Fagron	Neomycin sulfate	14A31-B02-294696	4
Fagron	Neomycin sulfate	4312I-01926	1
Fagron	Neomycin sulfate	352954	1
Fagron	Neomycin sulfate	16B28-B05-352957	1
Fagron	Neomycin sulfate	14A31B02294696	1
Fagron	Neomycin sulfate	14A31-B02-294697	1
Fagron	Neomycin sulfate	14A31B02	2
Fagron	Neomycin sulfate	14A31-B02	1
Fagron	Neomycin sulfate	14A31-B02-296724	1
Fagron	Neomycin sulfate	18G28-B05-352957	2
Fagron	Neomycin sulfate	14A31-B02-302519	1
Fagron	Neomycin sulfate	14c27-B05-302516	1
Fagron	Neomycin sulfate	16A20-B10	1
Fagron	Neomycin sulfate	18G28-B05-352953	2
Fagron	Neomycin sulfate	14A31-B02-302519	3
Fagron	Neomycin sulfate	F08469	1
Fagron	Neomycin sulfate	16B10-B50	3
Fagron	Neomycin sulfate	100810-112-91	1
Fagron	Neomycin sulfate	10121-N01	1
Fagron	Neomycin sulfate	17F12-B04	1
Fagron	Neomycin sulfate	17G05-B04	1
Fagron	Neomycin sulfate	17B14-B08-334000	1
Fagron	Neomycin sulfate	17B14-B08	1
Fagron	Neomycin sulfate	18A25-B05-346136	1
Fagron	Neomycin sulfate	12I19-N14	1
Fagron	Nicotinamide	L14020084	1
Fagron	Nicotinamide	12J02-B03	1
Fagron	Nitrofurantoin	15J15-B06-317250	2
Fagron	Nitrofurantoin	15J15-B06-319031	1
Fagron	Nitrofurantoin	13D10-N03	1
Fagron	Nitrofurantoin	3071403	1
Fagron	Nitrofurantoin	11F07-N06	3
Fagron	Nitrofurantoin	71/14	1
Fagron	Nitrofurantoin	201203080296	1
Fagron	Nitrofurantoin	15J15-B06	2
Fagron	Nitrofurantoin	12C30-N03	1
Fagron	Nitrofurantoin	13A08-N01	1
Fagron	Nitrofurantoin	15J15-B06-315225	2
Fagron	Nitrofurantoin	15J15-B06-313939	1
Fagron	Nitrofurantoin	12F19-N02	1
Fagron	Nitrofurantoin	2714I-03527	1
Fagron	Nitrofurantoin	13B28-N07	1
Fagron	Norethisterone acetate	14H29-B03-321322	1
Fagron	Norethisterone acetate	18C12-B03-354815	1
Fagron	Norethisterone acetate	18C12-B03-349111	1
Fagron	Norethisterone acetate	17B21-B01-341608	1
Fagron	Norethisterone acetate	17B21-B01-341607	2
Fagron	Norethisterone acetate	13D19-N06	1
Fagron	Norethisterone acetate	14D08-B01-292738	1
Fagron	Norethisterone acetate	17B21-B01	1
Fagron	Norethisterone acetate	14H29-B03-302277	2
Fagron	Norethisterone acetate	14h29-b03-310697	1
Fagron	Norethisterone acetate	14H29-B03-299313	1
Fagron	Norethisterone acetate	12C06-N10	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Norethisterone acetate	22061206	1
Fagron	Norethisterone acetate	14d08-b01-292736	1
Fagron	Norethisterone acetate	14H29-B03-325224	1
Fagron	Nystatin	4014376	1
Fagron	Nystatin	13C11-N08	1
Fagron	Nystatin	20121202	1
Fagron	Nystatin	12G04-N12	1
Fagron	Nystatin	11K11-No4	1
Fagron	Nystatin	12C29-N09	1
Fagron	Nystatin	11C21-N10	1
Fagron	Nystatin	17J25-B01-348043	1
Fagron	Octenidine dihydrochl...	187223403	1
Fagron	Ofloxacin	13A17-N01	2
Fagron	Ofloxacin	15H19-B01-314470	1
Fagron	Omeprazole	17K28-B07-348698	4
Fagron	Omeprazole	15B26-B40-305669	2
Fagron	Omeprazole	15K12-B02-314608	3
Fagron	Omeprazole	15K12-B02	1
Fagron	Omeprazole	15k12-B02	1
Fagron	Omeprazole	15K12-B02-525174	2
Fagron	Omeprazole	17B28-B02-525557	1
Fagron	Omeprazole	OPFP16030	1
Fagron	Omeprazole	18E22B03353496	1
Fagron	Omeprazole	16I20-B03-328940	1
Fagron	Omeprazole	16B16-B01-52535	1
Fagron	Omeprazole	17h15b0334821	1
Fagron	Omeprazole	17h15-b03-342821	1
Fagron	Omeprazole	16/20-B03-328940	3
Fagron	Omeprazole	17B28-B02-525511	1
Fagron	Omeprazole	17H15-B03-340200	1
Fagron	Omeprazole	16I20-B03-328940	1
Fagron	Ornithine hydrochloride	13A03-N05	1
Fagron	Ornithine hydrochloride	15A13-B03	1
Fagron	Orris root	11026705	3
Fagron	Oxybuprocaine hydroch...	13C01-N01	4
Fagron	Oxybuprocaine hydroch...	13J24B07295015	1
Fagron	Oxybuprocaine hydroch...	13h13-n05	2
Fagron	Oxybuprocaine hydroch...	4612M-03639	1
Fagron	Oxybuprocaine hydroch...	12L19N01	1
Fagron	Oxybuprocaine hydroch...	3613N-07321	1
Fagron	Oxybuprocaine hydroch...	15A26-B07	1
Fagron	Oxybutynin hydrochlor...	13H26-B03-287325	1
Fagron	Oxybutynin hydrochlor...	13H26-B03	1
Fagron	Oxybutynin hydrochlor...	14A29-B09-307322	5
Fagron	Oxybutynin hydrochlor...	14a29-b09-307322	2
Fagron	Oxybutynin hydrochlor...	127979	1
Fagron	Oxybutynin hydrochlor...	14/16_B03_300207	1
Fagron	Oxybutynin hydrochlor...	14A29_B09_307322	2
Fagron	Oxybutynin hydrochlor...	15F05-B03-315221	5
Fagron	Oxybutynin hydrochlor...	15F05-b03-315220	1
Fagron	Oxybutynin hydrochlor...	14I16-B03-300207	1
Fagron	Oxybutynin hydrochlor...	14/16-B03-300207	4
Fagron	Oxybutynin hydrochlor...	1709250/16I23-B03-33	1
Fagron	Oxybutynin hydrochlor...	16I23-b03-328916	1
Fagron	Oxybutynin hydrochlor...	16/23-b03-333371	1
Fagron	Oxybutynin hydrochlor...	18A12-B03-346964	2
Fagron	Oxybutynin hydrochlor...	18A12-B03-344347	1
Fagron	Oxybutynin hydrochlor...	16I23-B03-328916	2
Fagron	Oxybutynin hydrochlor...	18A12-B03-350471	1
Fagron	Oxybutynin hydrochlor...	14I16-B03	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Oxybutynin hydrochlor...	16A18-B02-322538	3
Fagron	Oxybutynin hydrochlor...	14/16-B03-300209	1
Fagron	Oxybutynin hydrochlor...	1607080	1
Fagron	Oxybutynin hydrochlor...	16A18-B02-326183	4
Fagron	Oxybutynin hydrochlor...	18A12-B03-344348	1
Fagron	Oxytetracycline hydro...	14I22-B01-301464	1
Fagron	Oxytetracycline hydro...	13L19-B33-290731	3
Fagron	Oxytetracycline hydro...	912A-02865	1
Fagron	Oxytetracycline hydro...	13D08-N11	6
Fagron	Oxytetracycline hydro...	15L16-B02-322908	1
Fagron	Oxytetracycline hydro...	14/22-B01-301465	1
Fagron	Oxytetracycline hydro...	14I22-B01-301463	1
Fagron	Oxytetracycline hydro...	14I22-B01	1
Fagron	Oxytetracycline hydro...	13I26-N11	1
Fagron	Oxytetracycline hydro...	4-49732A-2	1
Fagron	Oxytetracycline hydro...	130602203	1
Fagron	Oxytetracycline hydro...	12A26-N05	4
Fagron	Oxytetracycline hydro...	C22H25CIN209	1
Fagron	Oxytetracycline hydro...	14I22-B01-301465	1
Fagron	Oxytetracycline hydro...	12G10-N04	8
Fagron	Oxytetracycline hydro...	14/22-B01-301464	1
Fagron	Oxytetracycline hydro...	14I22-B01-301463	1
Fagron	Oxytetracycline hydro...	15L16-b02-322907	1
Fagron	Oxytetracycline hydro...	15L16-B02-322907	1
Fagron	Oxytetracycline hydro...	15D29-B99-308608	1
Fagron	Oxytetracycline hydro...	13L19-B33	5
Fagron	Paracetamol	13L16-B02-322725	1
Fagron	Paraffin, hard	14L01-B06	1
Fagron	Peppermint leaf, cut	13030803	3
Fagron	Peppermint leaf, cut	1234567	2
Fagron	Pheniramine maleate	14C25-B01-307150	3
Fagron	Pheniramine maleate	14C25-B01-328732	2
Fagron	Pheniramine maleate	12I06-N02	1
Fagron	Pheniramine maleate	14c25b01307150	1
Fagron	Pheniramine maleate	14c25-b01-307150	1
Fagron	Pheniramine maleate	12L04-N05	4
Fagron	Pheniramine maleate	14C25-B01-292540	2
Fagron	Pheniramine maleate	14C25-B01-295293	3
Fagron	Phenylbutazon	11H26-N02	1
Fagron	Phenylephrine hydroch...	12J19-B40-288253	2
Fagron	Phenylephrine hydroch...	16A20-B07-323398	3
Fagron	Phenylephrine hydroch...	17C23-B04-334819	3
Fagron	Phenylephrine hydroch...	16A20-B07-317987	2
Fagron	Phenylephrine hydroch...	16A-B07-317987	1
Fagron	Phenylephrine hydroch...	12031801	1
Fagron	Phenylephrine hydroch...	17c23-b04-333880	1
Fagron	Phenylephrine hydroch...	15F09-B08-316796	1
Fagron	Phenylephrine hydroch...	12C08-N17	1
Fagron	Phenylephrine hydroch...	13A03-N09	3
Fagron	Phenylephrine hydroch...	12A03-N09	1
Fagron	Phenylephrine hydroch...	13K25-B02-293844	3
Fagron	Phenylephrine hydroch...	17H18-B02-345401	2
Fagron	Phenylephrine hydroch...	14H13-B02-299814	1
Fagron	Phenylephrine hydroch...	14H13-B02-299816	2
Fagron	Phenylephrine hydroch...	12g19-n04	1
Fagron	Phenylephrine hydroch...	12G19-N04	1
Fagron	Phenylephrine hydroch...	13B14-N03	1
Fagron	Phenylephrine hydroch...	10091504	1
Fagron	Phenylephrine hydroch...	13G10-N04	1
Fagron	Phenylephrine hydroch...	14H13B02299817	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Phenylephrine hydroch...	1021402	1
Fagron	Phenylephrine hydroch...	14H13B02-299817	1
Fagron	Phenylephrine hydroch...	13K25-B02-293845	1
Fagron	Phenylephrine hydroch...	14H13-B02-309892	3
Fagron	Phenylephrine hydroch...	17C23-B04-333878	3
Fagron	Phenylephrine hydroch...	17C23-B04-333880	1
Fagron	Phenytoin	12I07-N02	1
Fagron	Phenytoin	13G26-N02	1
Fagron	Phenytoin	12C27N13	1
Fagron	Pilocarpine hydrochlo...	13H24-N03	1
Fagron	Pilocarpine hydrochlo...	15B27-B03-309493	1
Fagron	Poloxamers 407	13G29-N01	1
Fagron	Poloxamers 407	WPWI610C	1
Fagron	Poloxamers 407	12K23-N03	1
Fagron	Poloxamers 407	13G29-N01	1
Fagron	Poloxamers 407	5315I-07027	1
Fagron	Poloxamers 407	315E-07027	1
Fagron	Poloxamers 407	15A26-B03-303503	1
Fagron	Poloxamers 407	5211Q-07027	1
Fagron	Poloxamers 407	914I-07027	1
Fagron	Polymyxin B sulfate	13G05-N06	1
Fagron	Polymyxin B sulfate	12D05-N01	1
Fagron	Polymyxin B sulfate	13G05-N06	1
Fagron	Polymyxin B sulfate	13/13-N05	1
Fagron	Polymyxin B sulfate	12G13-N05	1
Fagron	Polymyxin B sulfate	13H28-B02-287954	1
Fagron	Polymyxin B sulfate	14F04-B03-296721	1
Fagron	Polymyxin B sulfate	912A-03315	1
Fagron	Polymyxin B sulfate	15H14-B01-315692	1
Fagron	Polymyxin B sulfate	16C01-B04-322012	3
Fagron	Polymyxin B sulfate	14H12-B06-304727	1
Fagron	Polymyxin B sulfate	16C01-B04	1
Fagron	Polymyxin B sulfate	17H23-B01-344840	1
Fagron	Polymyxin B sulfate	16G21-B03-333808	1
Fagron	Potassium citrate	1234567	1
Fagron	Potassium hydrocarbon...	12K27-N06	1
Fagron	Potassium hydrocarbon...	120801001	1
Fagron	Potassium hydrocarbon...	12141911	1
Fagron	Potassium hydrocarbon...	13031205	1
Fagron	Potassium sorbate	3912A-01748	1
Fagron	Potassium sorbate	4314U-01748	1
Fagron	Potassium sorbate	3115M-01748	1
Fagron	Povidone 25, soluble	24701988Q0	1
Fagron	Povidone 25, soluble	15C12-B02-327866	1
Fagron	Povidone 25, soluble	11276509T0	1
Fagron	Prasterone	13J10-B02	1
Fagron	Prasterone	14C18-Bo2-293603	1
Fagron	Prasterone	14C18-B02-293602	1
Fagron	Prasterone	7013260004	1
Fagron	Prasterone	14c18b02293602	1
Fagron	Prasterone	14h12b05-297723	1
Fagron	Prasterone	13J10-B02-287177	1
Fagron	Prasterone	12K30-N06	3
Fagron	Prasterone	19061304	1
Fagron	Prasterone	12C26-N09	1
Fagron	Prasterone	14C18-B02	1
Fagron	Prasterone	14C18-B02-293606	1
Fagron	Prasterone	27071705	1
Fagron	Prasterone	DHEAA150401-1	1
Fagron	Prasterone	C04-140804(M)	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Prasterone	15G14-B07-310827	1
Fagron	Prasterone	15A15-B40-307323	2
Fagron	Prasterone	18C12	1
Fagron	Prasterone	18B02-B02-345147	1
Fagron	Prasterone	1421-B03-300211	1
Fagron	Prasterone	16B19-B04-322246	2
Fagron	Prasterone	17D11-B10-336395	1
Fagron	Prasterone	15I18-B09-316798	1
Fagron	Prasterone	16L13-B01-332104	1
Fagron	Prasterone	16J26-B03-327602	1
Fagron	Prasterone	16I02-B02-326485	2
Fagron	Prasterone	14H12-B05-297727	2
Fagron	Prasterone	18B13-B09_346424	1
Fagron	Prasterone	12G04-N19	3
Fagron	Prednicarbate, micron...	12D27-N04	4
Fagron	Prednicarbate, micron...	15g10-b02-310848	1
Fagron	Prednicarbate, micron...	15G10-B02-314082	6
Fagron	Prednicarbate, micron...	15G10-B02-310848	5
Fagron	Prednicarbate, micron...	NU0823M	1
Fagron	Prednicarbate, micron...	15D15-B03-30816	1
Fagron	Prednicarbate, micron...	15G10-B02	1
Fagron	Prednicarbate, micron...	15G10-B02-311819	1
Fagron	Prednicarbate, micron...	15J22-B04-316799	5
Fagron	Prednicarbate, micron...	15J22-B04-316799	1
Fagron	Prednicarbate, micron...	15D15-B03-309682	1
Fagron	Prednicarbate, micron...	15J22-B04-318528	5
Fagron	Prednicarbate, micron...	15J22-B04-318930	3
Fagron	Prednicarbate, micron...	15J-B04-316799	1
Fagron	Prednicarbate, micron...	160518P	1
Fagron	Prednicarbate, micron...	160314PC	1
Fagron	Prednicarbate, micron...	16A18/B01	1
Fagron	Prednicarbate, micron...	16E12-B01	3
Fagron	Prednicarbate, micron...	16E12-B01-325758	4
Fagron	Prednicarbate, micron...	16E12-B01-328270	3
Fagron	Prednicarbate, micron...	16E12-B01-327344	3
Fagron	Prednicarbate, micron...	NT0555M13I23-N04	1
Fagron	Prednicarbate, micron...	15D15-B03-308162	1
Fagron	Prednicarbate, micron...	15D15-B03-308161	6
Fagron	Prednicarbate, micron...	14A27-B02	1
Fagron	Prednicarbate, micron...	NU0823M	1
Fagron	Prednicarbate, micron...	13I23-N04	6
Fagron	Prednicarbate, micron...	13/23-N04	4
Fagron	Prednicarbate, micron...	13G05-N03	4
Fagron	Prednicarbate, micron...	9E04-N02	1
Fagron	Prednicarbate, micron...	9G14-N06	1
Fagron	Prednicarbate, micron...	13c25-no6	1
Fagron	Prednicarbate, micron...	13F11-N01	3
Fagron	Prednicarbate, micron...	7061303	1
Fagron	Prednicarbate, micron...	13C25-N02	1
Fagron	Prednicarbate, micron...	14A27-B02-291258	2
Fagron	Prednicarbate, micron...	16J27-b01-330976	1
Fagron	Prednicarbate, micron...	13372408	2
Fagron	Prednicarbate, micron...	12B20-N06	2
Fagron	Prednicarbate, micron...	14a27-b02-291256	2
Fagron	Prednicarbate, micron...	14B13-B31-291255	1
Fagron	Prednicarbate, micron...	11L23-N07	5
Fagron	Prednicarbate, micron...	14A27-B02-295796	7
Fagron	Prednicarbate, micron...	14A27-B02-295797	6
Fagron	Prednicarbate, micron...	14A27-B02-295793	1
Fagron	Prednicarbate, micron...	701867-0001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Prednicarbate, micron...	14J16B06303466	1
Fagron	Prednicarbate, micron...	14J16-B06-303465	10
Fagron	Prednicarbate, micron...	13C25-N06	6
Fagron	Prednicarbate, micron...	16E12-B01-330093	2
Fagron	Prednicarbate, micron...	14A27-B02-291256	14
Fagron	Prednicarbate, micron...	16J27-B01	2
Fagron	Prednicarbate, micron...	17J04-B05-343934	1
Fagron	Prednicarbate, micron...	NX1422M	1
Fagron	Prednicarbate, micron...	174J04-B05-343935	1
Fagron	Prednicarbate, micron...	17J04-B05-341872	7
Fagron	Prednicarbate, micron...	17c27-b01-339079	1
Fagron	Prednicarbate, micron...	16/27-B01-332958	2
Fagron	Prednicarbate, micron...	17J04-B05-346113	8
Fagron	Prednicarbate, micron...	17J04-B05-340570	2
Fagron	Prednicarbate, micron...	170706PC,NW1915M	1
Fagron	Prednicarbate, micron...	17C27-B01-335531	1
Fagron	Prednicarbate, micron...	16J27-B01-332958	8
Fagron	Prednicarbate, micron...	17C27-B01-335530	9
Fagron	Prednicarbate, micron...	16J27-B01-330977	1
Fagron	Prednicarbate, micron...	16j27b01330976	1
Fagron	Prednicarbate, micron...	17C27-B01-339079	4
Fagron	Prednicarbate, micron...	17J04-B05-326113	1
Fagron	Prednicarbate, micron...	18B27-B01-347486	9
Fagron	Prednicarbate, micron...	18B27-B01-0001	1
Fagron	Prednicarbate, micron...	19E14-B01-191431	1
Fagron	Prednicarbate, micron...	18E14-B03-355727	1
Fagron	Prednicarbate, micron...	18B27-B01-346748	1
Fagron	Prednicarbate, micron...	18E14-B03	1
Fagron	Prednicarbate, micron...	18E14-B03-355726	3
Fagron	Prednicarbate, micron...	160902201	2
Fagron	Prednicarbate, micron...	16J28-B01-330976	1
Fagron	Prednicarbate, micron...	18E14B03-352560	1
Fagron	Prednicarbate, micron...	18B27-B01-351177	1
Fagron	Prednicarbate, micron...	18E14-B03-352631	4
Fagron	Prednicarbate, micron...	18B27-01-350264	1
Fagron	Prednicarbate, micron...	1612-B01-325758	1
Fagron	Prednicarbate, micron...	18B27-B01_347486	1
Fagron	Prednicarbate, micron...	18B27-01-347486	1
Fagron	Prednicarbate, micron...	18B27-B01-350264	4
Fagron	Prednicarbate, micron...	16j27-b01-330976	1
Fagron	Prednisolon dihydroge...	11K11-N11	2
Fagron	Prednisolon dihydroge...	15K27-B09-318738	2
Fagron	Prednisolon dihydroge...	17I27-B03	2
Fagron	Prednisolon dihydroge...	15J01-B02-317180	2
Fagron	Prednisolon dihydroge...	15J01-B02	4
Fagron	Prednisolon dihydroge...	15J01-B02-312722	5
Fagron	Prednisolon dihydroge...	14g30-b09-303014	1
Fagron	Prednisolon dihydroge...	15J01-B02-313610	4
Fagron	Prednisolon dihydroge...	15f29-b11-309694	1
Fagron	Prednisolon dihydroge...	15K27-B09	3
Fagron	Prednisolon dihydroge...	15K27-B09-322163	3
Fagron	Prednisolon dihydroge...	NU1356	1
Fagron	Prednisolon dihydroge...	15K27-B09-326263	3
Fagron	Prednisolon dihydroge...	15K27B09-326263	1
Fagron	Prednisolon dihydroge...	16I14-B03-328918	2
Fagron	Prednisolon dihydroge...	14G30-B09	1
Fagron	Prednisolon dihydroge...	14G30-B09-303014	4
Fagron	Prednisolon dihydroge...	14G30-B09-303012	6
Fagron	Prednisolon dihydroge...	16/14-B03-328169	1
Fagron	Prednisolon dihydroge...	16I14-B03	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Prednisolon dihydroge...	14E07-B06-295251	3
Fagron	Prednisolon dihydroge...	17B08-B05-337466	3
Fagron	Prednisolon dihydroge...	13/27-N01	1
Fagron	Prednisolon dihydroge...	13I27-N01	3
Fagron	Prednisolon dihydroge...	13H19N04	1
Fagron	Prednisolon dihydroge...	13C26N06	1
Fagron	Prednisolon dihydroge...	15k27-b09-318738	1
Fagron	Prednisolon dihydroge...	17B08B05337466	1
Fagron	Prednisolon dihydroge...	17/27-b03-342560	1
Fagron	Prednisolon dihydroge...	15K27-B09-322162	1
Fagron	Prednisolon dihydroge...	18D23-B08-353343	1
Fagron	Prednisolon dihydroge...	15F29-B11-309694	1
Fagron	Prednisolon dihydroge...	14E07-B06-295253	1
Fagron	Prednisolon dihydroge...	12H30-N08	2
Fagron	Prednisolon dihydroge...	13A10-N04	3
Fagron	Prednisolon dihydroge...	11C18-N11	1
Fagron	Prednisolone 1:10 (ri...	12L19-N06	1
Fagron	Prednisolone 1:10 (ri...	12D04-N01	1
Fagron	Prednisolone acetate	200540-0013	1
Fagron	Prednisolone acetate	13K20-B04-293970	1
Fagron	Prednisolone acetate	15b03-b05-302564	1
Fagron	Prednisolone acetate	115B03-B05-308582	1
Fagron	Prednisolone acetate	14D04-B03-299294	5
Fagron	Prednisolone acetate	15F29-B07-321001	2
Fagron	Prednisolone acetate	11I09N-08	1
Fagron	Prednisolone acetate	12k20-b04293971	1
Fagron	Prednisolone acetate	11I09-N08	1
Fagron	Prednisolone acetate	12G25-N05	5
Fagron	Prednisolone acetate	11A31-N03	1
Fagron	Prednisolone acetate	13K20-B04-289494	2
Fagron	Prednisolone acetate	13K20-B04-293971	2
Fagron	Prednisolone acetate	14G30-B13	2
Fagron	Prednisolone acetate	12J10-N05	1
Fagron	Prednisolone acetate	15I23-B04-340335	1
Fagron	Prednisolone acetate	12H14-N07	5
Fagron	Prednisolone acetate	12L11-B04-286905	1
Fagron	Prednisolone acetate	15B03-B05-302566	1
Fagron	Prednisolone acetate	K07A20130903M	1
Fagron	Prednisolone acetate	15803-B05	1
Fagron	Prednisolone acetate	14D04B03299295	1
Fagron	Prednisolone acetate	15I23-B04-336299	1
Fagron	Prednisolone acetate	15I23-B04-341793	1
Fagron	Prednisolone acetate	12C19-N03	4
Fagron	Prednisolone acetate	6082013D	1
Fagron	Prednisolone acetate	12C19-N03	1
Fagron	Prednisolone acetate	14G30-B13-299293	2
Fagron	Prednisolone acetate	15B03-B05-302564	13
Fagron	Prednisolone acetate	14G30-B13-302562	1
Fagron	Prednisolone acetate	17189601	1
Fagron	Prednisolone acetate	12L11-B40	1
Fagron	Prednisolone acetate	15B03-B05-308582	3
Fagron	Prednisolone acetate	13K20-B04-289495	3
Fagron	Prednisolone acetate	13D03-N06	1
Fagron	Prednisolone acetate	15I23-B04-336299	1
Fagron	Prednisolone acetate	1523-B04	1
Fagron	Prednisolone acetate	5444	1
Fagron	Prednisolone acetate	15/23-B04-334688	1
Fagron	Prednisolone acetate	15B03-B05	1
Fagron	Prednisolone acetate	15F29-B07-323268	1
Fagron	Prednisolone acetate	15I23-B04-327929	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Prednisolone acetate	15I23-B04-331154	3
Fagron	Prednisolone acetate	15F29-B07-317179	1
Fagron	Prednisolone acetate	15B03B05302564	1
Fagron	Prednisolone acetate	15/23-B04-341793	1
Fagron	Prednisolone acetate	15I23-B4-334688	1
Fagron	Prednisolone acetate	15I23-B04-334702	1
Fagron	Prednisolone acetate	12h14-n07	1
Fagron	Prednisolone acetate	113Q-02106	1
Fagron	Prednisolone acetate	15/23-B04-334702	1
Fagron	Prednisolone acetate	13A07-N07	9
Fagron	Prednisolone acetate	12J10/N05	1
Fagron	Prednisolone acetate	180713a	1
Fagron	Prednisolone acetate	12L11-B40-286905	14
Fagron	Prednisolone acetate	12L11-B40-286915	2
Fagron	Prednisolone acetate	12L05-N01	6
Fagron	Prednisolone acetate	13a07n07	1
Fagron	Prednisolone hydrogen...	14G29-B07-304733	1
Fagron	Prednisolone hydrogen...	15C11-B05-344778	3
Fagron	Prednisolone hydrogen...	15C11-B05-338297	1
Fagron	Prednisolone hydrogen...	12L14-N04	1
Fagron	Prednisolone hydrogen...	15c11-B05338297	1
Fagron	Prednisolone, microni...	13J28-B06-289670	3
Fagron	Prednisolone, microni...	613Q-02105	1
Fagron	Prednisolone, microni...	17G2B03-345538	1
Fagron	Prednisolone, microni...	17g2828-b03-345539	1
Fagron	Prednisolone, microni...	15D14-B01	1
Fagron	Prednisolone, microni...	17D14-Bo1-338573	1
Fagron	Prednisolone, microni...	17G28-B03-345538	2
Fagron	Prednisolone, microni...	17G28-B03-345539	1
Fagron	Prednisolone, microni...	15D14-B01-308425	3
Fagron	Prednisolone, microni...	14/10-B03-302271	1
Fagron	Prednisolone, microni...	14 10-B03-302265	1
Fagron	Prednisolone, microni...	13H07-B02-285590	1
Fagron	Prednisolone, microni...	18C19-B07-351226	2
Fagron	Prednisolone, microni...	14/10B03302265	1
Fagron	Prednisolone, microni...	17014-801-340448	1
Fagron	Prednisolone, microni...	14B18-B05294966	1
Fagron	Prednisolone, microni...	14B18-B05-294969	2
Fagron	Prednisolone, microni...	13H07-B02-285592	1
Fagron	Prednisolone, microni...	14I10-B03302265	1
Fagron	Prednisolone, microni...	13H07-B02	1
Fagron	Prednisolone, microni...	1410-B03	1
Fagron	Prednisolone, microni...	14B18-B05-294966	1
Fagron	Prednisolone, microni...	18H07-B04-355631	1
Fagron	Prednisolone, microni...	K04C20130203M	1
Fagron	Prednisolone, microni...	14L10-B03-302269	1
Fagron	Prednisolone, microni...	14B18-Bo5-294966	1
Fagron	Prednisolone, microni...	14I10-B03-302265	2
Fagron	Prednisolone, microni...	18A29-B02-348540	1
Fagron	Prednisolone, microni...	18H07-B04-353326	2
Fagron	Prednisolone, microni...	130F-N04	1
Fagron	Prednisolone, microni...	18h07-b04-353326	1
Fagron	Prednisolone, microni...	18H07-B04	2
Fagron	Prednisolone, microni...	18C19-B07	1
Fagron	Prednisolone, microni...	K04C20140502	1
Fagron	Prednisolone, microni...	18H07-804	1
Fagron	Prednisolone, microni...	14I10-B03	1
Fagron	Prednisolone, microni...	18H07-B04-354523	2
Fagron	Prednisolone, microni...	13H07-B02	1
Fagron	Prednisolone, microni...	K04C20131102	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Prednisolone, microni...	13H7-B02-285590	1
Fagron	Prednisolone, microni...	13H07-B02-285593	2
Fagron	Prednisolone, microni...	18A29-B02	1
Fagron	Prednisolone, microni...	18C19-B07-353242	1
Fagron	Prednisolone, microni...	14B18-B05	1
Fagron	Prednisolone, microni...	16C15-B01-324942	6
Fagron	Prednisolone, microni...	16c15-B01-320761	1
Fagron	Prednisolone, microni...	1711021A	1
Fagron	Prednisolone, microni...	16E09-B075411	1
Fagron	Prednisolone, microni...	K04C20121006	1
Fagron	Prednisolone, microni...	13E13-N06	1
Fagron	Prednisolone, microni...	15B18-B05-294969	1
Fagron	Prednisolone, microni...	13H07-B02-285590	1
Fagron	Prednisolone, microni...	13H07-B020	1
Fagron	Prednisolone, microni...	13E15-N06	4
Fagron	Prednisolone, microni...	16C15-B01-321979	1
Fagron	Prednisolone, microni...	14/10-B03-302265	3
Fagron	Prednisolone, microni...	13H07-B02-285593	1
Fagron	Prednisolone, microni...	13DE15-N06	1
Fagron	Prednisolone, microni...	14I10-B03-302269	2
Fagron	Prednisolone, microni...	K04C20121006	1
Fagron	Prednisolone, microni...	16A20-B03-319015	7
Fagron	Prednisolone, microni...	13E15-N06	1
Fagron	Prednisolone, microni...	4313E-02105	1
Fagron	Prednisolone, microni...	NEPL170210	1
Fagron	Prednisolone, microni...	16C15-B02-324945	1
Fagron	Prednisolone, microni...	16c15-b02	1
Fagron	Prednisolone, microni...	16E09/B07/328220	1
Fagron	Prednisolone, microni...	16E09-B07_328220	1
Fagron	Prednisolone, microni...	13E15-N06	1
Fagron	Prednisolone, microni...	1717M-02105	1
Fagron	Prednisolone, microni...	13A07N14	1
Fagron	Prednisolone, microni...	16E09-B07-328220	7
Fagron	Prednisolone, microni...	16E09-B07-328150	12
Fagron	Prednisolone, microni...	12K14-N03	8
Fagron	Prednisolone, microni...	16A20-B03	4
Fagron	Prednisolone, microni...	16H04-B03-332929	9
Fagron	Prednisolone, microni...	12c22-N05	1
Fagron	Prednisolone, microni...	15K03B04314226	1
Fagron	Prednisolone, microni...	15K03_B04_314226	1
Fagron	Prednisolone, microni...	16C15-B02	1
Fagron	Prednisolone, microni...	16c15-802	1
Fagron	Prednisolone, microni...	17D14-B01-338573	2
Fagron	Prednisolone, microni...	13H07-B02-285592	2
Fagron	Prednisolone, microni...	13J28-B06	2
Fagron	Prednisolone, microni...	16A20-B03-316421	2
Fagron	Prednisolone, microni...	16e09-B07-329686	1
Fagron	Prednisolone, microni...	15K03-B04-314226	13
Fagron	Prednisolone, microni...	15K03-B04-314358	7
Fagron	Prednisolone, microni...	1410-B03-302269	2
Fagron	Prednisolone, microni...	12D30-N04	11
Fagron	Prednisolone, microni...	15K03-B04-314357	2
Fagron	Prednisolone, microni...	17F07-B05-341575	3
Fagron	Prednisolone, microni...	14/10-B03	1
Fagron	Prednisolone, microni...	150528P	1
Fagron	Prednisolone, microni...	14I10-B03-302266	1
Fagron	Prednisolone, microni...	15K03-B04-314226	1
Fagron	Prednisolone, microni...	16C15-B01-321979	13
Fagron	Prednisolone, microni...	11J26-N04	2
Fagron	Prednisolone, microni...	17F07-B05	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Prednisolone, microni...	11D13-N05	1
Fagron	Prednisolone, microni...	16A20-B03-316420	6
Fagron	Prednisolone, microni...	16C15B01320761	1
Fagron	Prednisolone, microni...	13A07-N14	1
Fagron	Prednisolone, microni...	16C15-B01-319992	2
Fagron	Prednisolone, microni...	11G22-N08	1
Fagron	Prednisolone, microni...	16C15B01321979	1
Fagron	Prednisolone, microni...	17D14-B01	1
Fagron	Prednisolone, microni...	15k03-b04-314226	1
Fagron	Prednisolone, microni...	16C15-B02-324943	6
Fagron	Prednisolone, microni...	12D30-N04	1
Fagron	Prednisolone, microni...	NEPL151202	1
Fagron	Prednisolone, microni...	13H07-B02-285593	1
Fagron	Prednisolone, microni...	NEPL170209	1
Fagron	Prednisone	13C13-N08	2
Fagron	Prednisone	15218608	1
Fagron	Prednisone	12E02-N01	1
Fagron	Prednisone	12E02-N10	1
Fagron	Prednisone	21021401	1
Fagron	Prednisone	16H04-B04-334368	1
Fagron	Pregnenolone	16L29-B02331611	1
Fagron	Pregnenolone	14B18-B04-291917	1
Fagron	Prilocaine hydrochlor...	15L21-B02325756	1
Fagron	Prilocaine hydrochlor...	15L21-B02-332956	2
Fagron	Prilocaine hydrochlor...	17Ho2Bo5348769	1
Fagron	Prilocaine hydrochlor...	13I032-N03	1
Fagron	Prilocaine hydrochlor...	15E15-B04-307730	1
Fagron	Prilocaine hydrochlor...	17H02-B05-340539	1
Fagron	Prilocaine hydrochlor...	15L21-B02_332956	1
Fagron	Prilocaine hydrochlor...	15L21-B02332956	1
Fagron	Prilocaine hydrochlor...	14H13-B03-297591	1
Fagron	Procaine hydrochloride	12I10-N07	1
Fagron	Procaine hydrochloride	13F12-B40	1
Fagron	Procaine hydrochloride	10H11-N25	1
Fagron	Procaine hydrochloride	1651212	1
Fagron	Procaine hydrochloride	13I09-N08	1
Fagron	Procaine hydrochloride	13/30-N06	1
Fagron	Progesterone, microni...	12D05-N06	3
Fagron	Progesterone, microni...	5013I-03235	1
Fagron	Progesterone, microni...	14D22-B04	2
Fagron	Progesterone, microni...	14c10-B01-293478	1
Fagron	Progesterone, microni...	14D22-Bo4	1
Fagron	Progesterone, microni...	1910213235	1
Fagron	Progesterone, microni...	12G05-N17	6
Fagron	Progesterone, microni...	12K05-N01	2
Fagron	Progesterone, microni...	13K28-B30-288261	1
Fagron	Progesterone, microni...	13K28-B30-288262	1
Fagron	Progesterone, microni...	1311-N11	1
Fagron	Progesterone, microni...	13K-B30-288261	1
Fagron	Progesterone, microni...	13H01-N03	6
Fagron	Progesterone, microni...	13/11-N11	1
Fagron	Progesterone, microni...	13b08-n01	1
Fagron	Progesterone, microni...	13G11-N02	5
Fagron	Progesterone, microni...	13A02-N01	5
Fagron	Progesterone, microni...	12k05-n07	1
Fagron	Progesterone, microni...	14C10B0129347	1
Fagron	Progesterone, microni...	1-111218	1
Fagron	Progesterone, microni...	12K05-N07	1
Fagron	Progesterone, microni...	13A02N01	1
Fagron	Progesterone, microni...	12A19-N04	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Progesterone, microni...	14022-B04-293476	1
Fagron	Progesterone, microni...	C01-130301	1
Fagron	Progesterone, microni...	16G15-B04-325136	3
Fagron	Progesterone, microni...	16A20-B02-321374	4
Fagron	Progesterone, microni...	16A20-B02-319103,201	1
Fagron	Progesterone, microni...	16G15-B04-326487	3
Fagron	Progesterone, microni...	16G15-B04-32579	1
Fagron	Progesterone, microni...	16G15-B04-325135	1
Fagron	Progesterone, microni...	16g115-b04-325136	1
Fagron	Progesterone, microni...	16G15-B04325136	1
Fagron	Progesterone, microni...	20170427Ma6	1
Fagron	Progesterone, microni...	16G15B04-325279	1
Fagron	Progesterone, microni...	16g15-b04	1
Fagron	Progesterone, microni...	EU-HTT160705	2
Fagron	Progesterone, microni...	16J05-B04-328303	2
Fagron	Progesterone, microni...	16G15-B04	2
Fagron	Progesterone, microni...	16L29-B03	1
Fagron	Progesterone, microni...	375779	1
Fagron	Progesterone, microni...	200650-0004	2
Fagron	Progesterone, microni...	16111-B04-314184	1
Fagron	Progesterone, microni...	16L29-B03-332060	5
Fagron	Progesterone, microni...	16D13-B04-323636	1
Fagron	Progesterone, microni...	16G15-B04-325755	1
Fagron	Progesterone, microni...	16A20-B02-318667	3
Fagron	Progesterone, microni...	15C31-B06-316105	3
Fagron	Progesterone, microni...	200650-315224	1
Fagron	Progesterone, microni...	16A20-B02-321375	2
Fagron	Progesterone, microni...	1511-B05-314184	1
Fagron	Progesterone, microni...	15I11-B05	1
Fagron	Progesterone, microni...	16A20-B02-321373	2
Fagron	Progesterone, microni...	14K11-B04-300921	1
Fagron	Progesterone, microni...	1511-805	1
Fagron	Progesterone, microni...	16a20-b02-321375	1
Fagron	Progesterone, microni...	16a20-b02-318667	1
Fagron	Progesterone, microni...	1511B05-315224	1
Fagron	Progesterone, microni...	16G15-B04-325279	5
Fagron	Progesterone, microni...	160A20-B02-321373	1
Fagron	Progesterone, microni...	16G15-B04325279	1
Fagron	Progesterone, microni...	15I11-B05-313598	1
Fagron	Progesterone, microni...	17C15-B03-333672	2
Fagron	Progesterone, microni...	17C15-B03-333671	3
Fagron	Progesterone, microni...	17C15-B03-333671	1
Fagron	Progesterone, microni...	17I08-B04-344773	1
Fagron	Progesterone, microni...	18B06-B10-346701	1
Fagron	Progesterone, microni...	18B06-B10-346702	1
Fagron	Progesterone, microni...	C01-170814	1
Fagron	Progesterone, microni...	18B06-B10-347484	4
Fagron	Progesterone, microni...	18B06-B10-349775	2
Fagron	Progesterone, microni...	18806-B10-349153	1
Fagron	Progesterone, microni...	18F22-B04-353582	1
Fagron	Progesterone, microni...	18F22-B04353582-P201	1
Fagron	Progesterone, microni...	18F22-B04-351435	2
Fagron	Progesterone, microni...	18F22-B04-354900	2
Fagron	Progesterone, microni...	18F22-B04	1
Fagron	Progesterone, microni...	18F22-B04-352735	1
Fagron	Progesterone, microni...	18F22-B04-351145	1
Fagron	Progesterone, microni...	18F22-B04-355395	1
Fagron	Progesterone, microni...	ApC01-180128P	1
Fagron	Progesterone, microni...	17/21-B03-342135	2
Fagron	Progesterone, microni...	17/21-B03-342134	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Progesterone, microni...	intern20180619Ma3	1
Fagron	Progesterone, microni...	16129b03332060	1
Fagron	Progesterone, microni...	16L29-B03-302060	1
Fagron	Progesterone, microni...	17C15B03333671	2
Fagron	Progesterone, microni...	17D25-B06-339655	2
Fagron	Progesterone, microni...	16L29-B-332061	1
Fagron	Progesterone, microni...	17D25-B06-339651	1
Fagron	Progesterone, microni...	17/21-B03-341517	2
Fagron	Progesterone, microni...	HTT140503	1
Fagron	Progesterone, microni...	17/21-B03-341518	1
Fagron	Progesterone, microni...	17I21-803	1
Fagron	Progesterone, microni...	17I21-B03-342135	1
Fagron	Progesterone, microni...	17I21-B03-343284	2
Fagron	Progesterone, microni...	17I21-B03-341518	1
Fagron	Progesterone, microni...	1511-B05-314842	1
Fagron	Progesterone, microni...	17I21-B03	1
Fagron	Progesterone, microni...	17I08-B04	2
Fagron	Progesterone, microni...	16L29-B03-332505	4
Fagron	Progesterone, microni...	15C31-B06-313599	1
Fagron	Progesterone, microni...	14C10-B01-293478	3
Fagron	Progesterone, microni...	15B12-B04-304837	2
Fagron	Progesterone, microni...	15B12-Bo4-304838	1
Fagron	Progesterone, microni...	14K11B04300920	1
Fagron	Progesterone, microni...	HTT150201	3
Fagron	Progesterone, microni...	15C31-B06	1
Fagron	Progesterone, microni...	15C31-B06-309790	2
Fagron	Progesterone, microni...	13k28-b30-288261	1
Fagron	Progesterone, microni...	15I11-B05-314842	1
Fagron	Progesterone, microni...	14L15-B05-302292	2
Fagron	Progesterone, microni...	15C31-B06-308439	7
Fagron	Progesterone, microni...	14K11-B04	3
Fagron	Progesterone, microni...	14D22-B293479	1
Fagron	Progesterone, microni...	HTT140302	2
Fagron	Progesterone, microni...	14H13-B05-300216	7
Fagron	Progesterone, microni...	14K11-B04-300920	7
Fagron	Progesterone, microni...	14C10-B01	1
Fagron	Progesterone, microni...	14D22-B04-293477	1
Fagron	Progesterone, microni...	13h01-n03	1
Fagron	Progesterone, microni...	14c10-b01-293478	1
Fagron	Progesterone, microni...	HTT140207	1
Fagron	Progesterone, microni...	15C31-B06-308796	2
Fagron	Progesterone, microni...	14k11-b04	1
Fagron	Progesterone, microni...	15C31-Bo6-308796	1
Fagron	Progesterone, microni...	14D22B04293476, 20160	1
Fagron	Progesterone, microni...	14K11-B04_300921	1
Fagron	Progesterone, microni...	15C31B06308439	1
Fagron	Progesterone, microni...	HTT140413	1
Fagron	Propranolol hydrochlo...	4071401	1
Fagron	Propranolol hydrochlo...	15K11-B01-316401	14
Fagron	Propranolol hydrochlo...	402310015	1
Fagron	Propranolol hydrochlo...	12F29-N06	1
Fagron	Propranolol hydrochlo...	15F16B02312979	1
Fagron	Propranolol hydrochlo...	13E15-B40	2
Fagron	Propranolol hydrochlo...	3112U-03609	1
Fagron	Propranolol hydrochlo...	17K22-B03	1
Fagron	Propranolol hydrochlo...	M130104	1
Fagron	Propranolol hydrochlo...	13E15_B40-297021	1
Fagron	Propranolol hydrochlo...	12K06N07	1
Fagron	Propranolol hydrochlo...	28051508	1
Fagron	Propranolol hydrochlo...	13E15-B40-297021	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Propranolol hydrochlo...	11I19-N09	1
Fagron	Propranolol hydrochlo...	12K06-N07	8
Fagron	Propranolol hydrochlo...	13E15-B40-294297	3
Fagron	Propranolol hydrochlo...	13J13-N04	1
Fagron	Propranolol hydrochlo...	18D23-B09	3
Fagron	Propranolol hydrochlo...	13E15-B40	3
Fagron	Propranolol hydrochlo...	13E15-B40-289355	3
Fagron	Propranolol hydrochlo...	13E15-B40-294297	3
Fagron	Propranolol hydrochlo...	13I13-N04	1
Fagron	Propranolol hydrochlo...	5101	1
Fagron	Propranolol hydrochlo...	3050PRII	1
Fagron	Propranolol hydrochlo...	15K11-Bo1-316401	1
Fagron	Propranolol hydrochlo...	11L01-N05	2
Fagron	Propranolol hydrochlo...	17K-22-B03-345013	1
Fagron	Propranolol hydrochlo...	13E15-B40-289355	9
Fagron	Propranolol hydrochlo...	13I3-N04	2
Fagron	Propranolol hydrochlo...	9L09-N06	1
Fagron	Propranolol hydrochlo...	3050PRRII	2
Fagron	Propranolol hydrochlo...	13I13-N04	8
Fagron	Propranolol hydrochlo...	18041807	1
Fagron	Propranolol hydrochlo...	12F29-N06	4
Fagron	Propranolol hydrochlo...	17K22-B03-342731	7
Fagron	Propranolol hydrochlo...	M170602	1
Fagron	Propranolol hydrochlo...	12D27-N14	5
Fagron	Propranolol hydrochlo...	16E09-B09333442	1
Fagron	Propranolol hydrochlo...	15K11-B01-16401	1
Fagron	Propranolol hydrochlo...	17K22-B03-350048	4
Fagron	Propranolol hydrochlo...	700556-0001	1
Fagron	Propranolol hydrochlo...	16E09-B09-326637	16
Fagron	Propranolol hydrochlo...	16E09-B09-333442	8
Fagron	Propranolol hydrochlo...	10L20N08	2
Fagron	Propranolol hydrochlo...	15K11-B01	4
Fagron	Propranolol hydrochlo...	9CMK6AA71374	1
Fagron	Propranolol hydrochlo...	31101201	1
Fagron	Propranolol hydrochlo...	11D01-N11	1
Fagron	Propranolol hydrochlo...	18D23-B09-351774	10
Fagron	Propranolol hydrochlo...	12A27-N04	1
Fagron	Propranolol hydrochlo...	15K11-B01-320607	9
Fagron	Propranolol hydrochlo...	17K22-b03-350048	1
Fagron	Propranolol hydrochlo...	20250711-8	1
Fagron	Propranolol hydrochlo...	17K22-B03-30048	1
Fagron	Propranolol hydrochlo...	13 13-N04	1
Fagron	Propranolol hydrochlo...	3913M-03609	1
Fagron	Propranolol hydrochlo...	13J13-N04	2
Fagron	Propranolol hydrochlo...	15K11-B01-322717	7
Fagron	Propyl-4-hydroxybenzo...	12B13-N04	1
Fagron	Propyl-4-hydroxybenzo...	610M-01931	1
Fagron	Propyphenazone	17G03-B01-347588	2
Fagron	Propyphenazone	16E02-B08-327928	2
Fagron	Propyphenazone	13K20-B03-319137	3
Fagron	Propyphenazone	13K20B03319137	1
Fagron	Propyphenazone	13k20-b03-319137	1
Fagron	Propyphenazone	13K-B03-319137	1
Fagron	Propyphenazone	13K20-B03-311715	3
Fagron	Propyphenazone	13K20-B03	2
Fagron	Propyphenazone	13K20-B03-302291	1
Fagron	Propyphenazone	13B12-B02-297171	1
Fagron	Pyridoxine hydrochlor...	12e02-n02	1
Fagron	Pyridoxine hydrochlor...	L13080030	1
Fagron	Pyridoxine hydrochlor...	1608035-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Quinidine sulfate dihydrochloride	13G08-N02	1
Fagron	Quinine hydrochloride	10N129B	1
Fagron	Quinine hydrochloride	13G10-N07	1
Fagron	Quinoline yellow	100427	1
Fagron	Quinoline yellow	14D18-B02-295783	1
Fagron	Quinoline yellow	9878	1
Fagron	Quinoline yellow	12B08-N08	1
Fagron	Quinoline yellow	11F30-N04	1
Fagron	Quinoline yellow	12E24-N05	1
Fagron	Quinoline yellow	12L11-N07	1
Fagron	Quinoline yellow	11K08-N01	2
Fagron	Resocrinol	15K24-B02-317584	4
Fagron	Resocrinol	11/13-N02	3
Fagron	Resocrinol	1112013A	1
Fagron	Resocrinol	13E24-N10	10
Fagron	Resocrinol	17H03-B11-342399	1
Fagron	Resocrinol	10B01-N06	2
Fagron	Resocrinol	13I18-B09	2
Fagron	Resocrinol	13/18-B09-291534	1
Fagron	Resocrinol	14/19050	1
Fagron	Resocrinol	13I18-B09-291533	1
Fagron	Resocrinol	16D29-B01-332738	2
Fagron	Resocrinol	16d29b01332738	1
Fagron	Resocrinol	13I18-B09-291534	2
Fagron	Resocrinol	118b09291534	1
Fagron	Resocrinol	14G10-B06-296296	8
Fagron	Resocrinol	18B08-B05-347787	4
Fagron	Resocrinol	14G10-B06-296297	2
Fagron	Resocrinol	14/16062	1
Fagron	Resocrinol	12g24-N04	1
Fagron	Resocrinol	13/17071	1
Fagron	Resocrinol	2112I-02170	1
Fagron	Resocrinol	11/13-no2	1
Fagron	Resocrinol	16d29b01323046	1
Fagron	Resocrinol	12K30-N03	4
Fagron	Resocrinol	12G24-N004	1
Fagron	Resocrinol	11L19-N19	2
Fagron	Resocrinol	16D29-B01-329760	3
Fagron	Resocrinol	12G24-N04	13
Fagron	Resocrinol	16D29-B01-323046	1
Fagron	Resocrinol	11A19-N09	1
Fagron	Resocrinol	14J15-B02	1
Fagron	Resocrinol	14J15-B02-307749	5
Fagron	Resocrinol	14G10B06	1
Fagron	Resocrinol	15/18082	1
Fagron	Rhubarb, powdered	11043821	1
Fagron	Riboflavin	14B03B07292071	1
Fagron	Riboflavin	14K03-B04-304718	1
Fagron	Riboflavin	167720-0005	1
Fagron	Riboflavin	13J29-B04-287533	1
Fagron	Riboflavin	14E01-B10-296718	1
Fagron	Riboflavin	13J29-B04-287534	1
Fagron	Riboflavin	16B18_b06_319026	1
Fagron	Riboflavin	15I09-B02-315216	1
Fagron	Riboflavin	1701050-01	1
Fagron	Riboflavin	18J10-B01-356142	1
Fagron	Riboflavin	15G20-B13-310997	2
Fagron	Rice starch	14F18-B04-301616	1
Fagron	Rice starch	12L19-N06	1
Fagron	Rutoside trihydrate	L14050127	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Rutoside trihydrate	13K12-B06	1
Fagron	Rutoside trihydrate	13K12B06295483	1
Fagron	Rutoside trihydrate	13K12-B06295483	1
Fagron	Saccharin sodium	12A13-N03	2
Fagron	Saccharin sodium	13D19-N01	1
Fagron	Saccharin sodium	SH349/12	1
Fagron	Saccharin sodium	12K26-N09	1
Fagron	Saffron filaments	380A120416	1
Fagron	Salicylic acid	18E24-B03	1
Fagron	Salicylic acid	102455-0002	1
Fagron	Salicylic acid	13B14-N01	20
Fagron	Salicylic acid	xxxxxxx	1
Fagron	Salicylic acid	12H10-N03	6
Fagron	Salicylic acid	12H10-N08	8
Fagron	Salicylic acid	3111Q-01045	1
Fagron	Salicylic acid	11197121	1
Fagron	Salicylic acid	12C20-N12	6
Fagron	Salicylic acid	11K10-N01	2
Fagron	Salicylic acid	16I15-B08	2
Fagron	Salicylic acid	16I15-B08-330773	1
Fagron	Salicylic acid	16/15-B08-332098	3
Fagron	Salicylic acid	16I15-B08-333932	2
Fagron	Salicylic acid	16I15B08330773	1
Fagron	Salicylic acid	16A18B04322905	2
Fagron	Salicylic acid	16I15-B08-327353	1
Fagron	Salicylic acid	17C21-B09	2
Fagron	Salicylic acid	17C21-B09-337173	3
Fagron	Salicylic acid	17C21-B09-337172	2
Fagron	Salicylic acid	17C21B09	1
Fagron	Salicylic acid	16I15-B08-332098	1
Fagron	Salicylic acid	17C21-B09-337860	2
Fagron	Salicylic acid	17c21-b09-337860	1
Fagron	Salicylic acid	18F05-B02	2
Fagron	Salicylic acid	16/15-B08-333933	1
Fagron	Salicylic acid	18F05-B02-350587	1
Fagron	Salicylic acid	18F05-B02-351974	1
Fagron	Salicylic acid	18E24-B03-354311	1
Fagron	Salicylic acid	18F05B02350587	2
Fagron	Salicylic acid	14L02B02303558	1
Fagron	Salicylic acid	16/15-B08-330773	3
Fagron	Salicylic acid	13H01-N01	7
Fagron	Salicylic acid	14E26-B04_294623	1
Fagron	Salicylic acid	14e26-b04-294623	1
Fagron	Salicylic acid	14E26-B04-294622	3
Fagron	Salicylic acid	16A18-B04322904	1
Fagron	Salicylic acid	14L02-B02	4
Fagron	Salicylic acid	14B05B30	1
Fagron	Salicylic acid	14E26B04	1
Fagron	Salicylic acid	14E26_B04	1
Fagron	Salicylic acid	16A18-B04	1
Fagron	Salicylic acid	RAS1135100	2
Fagron	Salicylic acid	2201115	1
Fagron	Salicylic acid	16A18-B04-322905	2
Fagron	Salicylic acid	12c20n12	1
Fagron	Salicylic acid	1140616	1
Fagron	Salicylic acid	14L02-B02-303557	10
Fagron	Salicylic acid	11L19-N21	1
Fagron	Salicylic acid	14E26B04294623	1
Fagron	Salicylic acid	14e26-B04-294623	1
Fagron	Salicylic acid	13H01-N08	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Salicylic acid	13H01N01	1
Fagron	Salicylic acid	13H01-N04	16
Fagron	Salicylic acid	RAS1135300	1
Fagron	Salicylic acid	13H01N04	4
Fagron	Salicylic acid	13K08-B31-290852	5
Fagron	Salicylic acid	14B05-B30-290854	4
Fagron	Salicylic acid	13K08-B31-290853	5
Fagron	Salicylic acid	16 15-B08-327353	1
Fagron	Salicylic acid	16A18-B04-322904	2
Fagron	Salicylic acid	290852	1
Fagron	Salicylic acid	14B05-B30	2
Fagron	Salicylic acid	14E26-B04-294623	6
Fagron	Salicylic acid	14E26-B04-294624	6
Fagron	Salicylic acid	14E26-B04	6
Fagron	Salicylic acid	12C20-N04	1
Fagron	Salicylic acid	14L02B02303557	1
Fagron	Salicylic acid	14E26-B04-296424	1
Fagron	Senna leaf, cut	314Q1514	1
Fagron	Sildenafil citrate	16F30-B02-331450	3
Fagron	Sildenafil citrate	16A18-B08	2
Fagron	Sildenafil citrate	16049817	1
Fagron	Sildenafil citrate	16F30-B02-324715	1
Fagron	Sildenafil citrate	17H31-B02-342459	3
Fagron	Sildenafil citrate	14H11-B01-303183	1
Fagron	Sildenafil citrate	16F30-B02	1
Fagron	Sildenafil citrate	16F30-B02-334977	1
Fagron	Sildenafil citrate	15F11-B02-320186	1
Fagron	Sildenafil citrate	16F30-B02-331449	1
Fagron	Sildenafil citrate	13J15-N01	2
Fagron	Sildenafil citrate	17071505	1
Fagron	Sildenafil citrate	15f11-b02-320186	1
Fagron	Sildenafil citrate	16F30-B02-324715	1
Fagron	Sildenafil citrate	15F11-B02	2
Fagron	Sodium alginate	14G22-B99	1
Fagron	Sodium ascorbate	12L27-N03	1
Fagron	Sodium ascorbate	13C15-N01	1
Fagron	Sodium benzoate	13C14-N01	3
Fagron	Sodium benzoate	13I27-B03-288806	4
Fagron	Sodium benzoate	24011401	1
Fagron	Sodium benzoate	15B16-B07-308570	1
Fagron	Sodium benzoate	1412	1
Fagron	Sodium benzoate	15363403	1
Fagron	Sodium benzoate	16B23-B06-322528	1
Fagron	Sodium benzoate	12J23-N05	4
Fagron	Sodium benzoate	12C05-N04	1
Fagron	Sodium benzoate	13C14-N01	1
Fagron	Sodium benzoate	test	1
Fagron	Sodium benzoate	15B16-B07	1
Fagron	Sodium benzoate	13I27-B03288805	1
Fagron	Sodium benzoate	14E12B07296717	1
Fagron	Sodium benzoate	12H08-N06	1
Fagron	Sodium benzoate	13/27b03288805	1
Fagron	Sodium benzoate	13127-B03-288806	1
Fagron	Sodium benzoate	13/27B03288805	1
Fagron	Sodium benzoate	13/27-B03-288805	1
Fagron	Sodium benzoate	13C14-N01	1
Fagron	Sodium benzoate	13H01-N02	1
Fagron	Sodium benzoate	12G04-N04	1
Fagron	Sodium benzoate	12D02-N06	1
Fagron	Sodium benzoate	14E12B07296717	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Sodium benzoate	14E12-B07-306570	1
Fagron	Sodium benzoate	14E12-B07	1
Fagron	Sodium benzoate	15B16-B07-306571	2
Fagron	Sodium benzoate	14E12Bo7296717	1
Fagron	Sodium citrate	815Q-01900	1
Fagron	Sodium citrate	15C27-B01-316106	1
Fagron	Sodium citrate	15c27-b01-321097	1
Fagron	Sodium citrate	15C27-B01-321097	1
Fagron	Sodium citrate	15C27-B01-312231	5
Fagron	Sodium citrate	15C27-B01-3211097	1
Fagron	Sodium citrate	13A07-N01	1
Fagron	Sodium citrate	14A13-B02-300395	4
Fagron	Sodium citrate	1176767	3
Fagron	Sodium citrate	14A13-B02-294541	8
Fagron	Sodium citrate	1172108	1
Fagron	Sodium citrate	13I24-B13-287921	1
Fagron	Sodium citrate	13B28-N02	1
Fagron	Sodium citrate	13A07-N01	1
Fagron	Sodium citrate	160614	1
Fagron	Sodium citrate	13D09-N01	4
Fagron	Sodium citrate	24121301	1
Fagron	Sodium citrate	Probe	1
Fagron	Sodium citrate	13B28-N02	6
Fagron	Sodium citrate	1161207	4
Fagron	Sodium citrate	12339502	1
Fagron	Sodium citrate	12E04-N05	1
Fagron	Sodium citrate	13A07-N01	5
Fagron	Sodium citrate	712A-01900	1
Fagron	Sodium citrate	1234567	1
Fagron	Sodium citrate	12J26-N02	2
Fagron	Sodium citrate	3311E-01900	1
Fagron	Sodium citrate	12I10-N01	1
Fagron	Sodium citrate	12A05-N04	2
Fagron	Sodium citrate	31051202	1
Fagron	Sodium citrate	13D09N01	1
Fagron	Sodium citrate	13/24-B13-287921	1
Fagron	Sodium citrate	179915	2
Fagron	Sodium citrate	13I24-B13-287921	1
Fagron	Sodium citrate	15C27-B01-309893	5
Fagron	Sodium citrate	14A13_b02_300395	1
Fagron	Sodium citrate	14369402	2
Fagron	Sodium citrate	17A26-B04-331528	14
Fagron	Sodium citrate	18C05-B03	1
Fagron	Sodium citrate	16F10-B01-327510	6
Fagron	Sodium citrate	17H07	1
Fagron	Sodium citrate	18B14-B08-348067	1
Fagron	Sodium citrate	18B14-B08	1
Fagron	Sodium citrate	18G25-B07-352737	1
Fagron	Sodium citrate	18G25-B07-354812	1
Fagron	Sodium citrate	1219546	1
Fagron	Sodium citrate	18L13B023549449	1
Fagron	Sodium citrate	17H07-B04-341746	2
Fagron	Sodium citrate	1190802	1
Fagron	Sodium citrate	17A26-B04-335883	2
Fagron	Sodium citrate	16F10-B1-327510	1
Fagron	Sodium citrate	15C27	1
Fagron	Sodium citrate	17H07-B04_338530	28
Fagron	Sodium citrate	17A26-B04-335680	24
Fagron	Sodium citrate	15C27-B01-323451	1
Fagron	Sodium citrate	17A26-B04-332871	26

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Sodium citrate	1211142	2
Fagron	Sodium citrate	16F10-B01-326088	3
Fagron	Sodium citrate	18C05-B03-350564	1
Fagron	Sodium citrate	17A26-B04-331525	8
Fagron	Sodium cromoglicate	13E22-N11	1
Fagron	Sodium dihydrogen pho...	12L05-N05	1
Fagron	Sodium dodecyl sulfate	16229205	4
Fagron	Sodium edetate	17D25-B08	3
Fagron	Sodium edetate	13/23-N03	4
Fagron	Sodium edetate	17D25-B08-334827	4
Fagron	Sodium edetate	17D25-B08-338202	1
Fagron	Sodium edetate	12E16-N04	3
Fagron	Sodium edetate	601H9118	1
Fagron	Sodium edetate	17D25B08343253	1
Fagron	Sodium edetate	13I23-N03	8
Fagron	Sodium edetate	18H03-B04	2
Fagron	Sodium edetate	18H03-B04-353034	5
Fagron	Sodium edetate	13E15-N05	16
Fagron	Sodium edetate	12g25-n12	1
Fagron	Sodium edetate	601HA062	3
Fagron	Sodium edetate	12E12-N04	1
Fagron	Sodium edetate	17D25-B08-343253	1
Fagron	Sodium edetate	18H03-B04-355280	2
Fagron	Sodium edetate	813Q-01925	1
Fagron	Sodium edetate	18h03-b04-356574	1
Fagron	Sodium edetate	601H1300	3
Fagron	Sodium edetate	17D25-B08-334827	1
Fagron	Sodium edetate	14G03B05309971	1
Fagron	Sodium edetate	13E15-B10-289034	1
Fagron	Sodium edetate	13G15-N03	4
Fagron	Sodium edetate	14G03-B05_307737	1
Fagron	Sodium edetate	13I23-N03	2
Fagron	Sodium edetate	601EB174	1
Fagron	Sodium edetate	14G03-B05-307737	1
Fagron	Sodium edetate	14G03-B05-307737	7
Fagron	Sodium edetate	14G03-B05-309971	1
Fagron	Sodium edetate	13K19B02290051	1
Fagron	Sodium edetate	14G03-B05-303420	1
Fagron	Sodium edetate	13K19-B02-290052	3
Fagron	Sodium edetate	13K19-B02	1
Fagron	Sodium edetate	14G03-B05	4
Fagron	Sodium edetate	13K19B02	1
Fagron	Sodium edetate	14G03-B05-30773	1
Fagron	Sodium edetate	11A31-N06Retest	1
Fagron	Sodium edetate	215I-01925	1
Fagron	Sodium edetate	601EB262	1
Fagron	Sodium edetate	601HA062	3
Fagron	Sodium edetate	5315A-01925	1
Fagron	Sodium edetate	16F21-B02-325891	2
Fagron	Sodium edetate	15D01-B05-317369	1
Fagron	Sodium edetate	601EE047	2
Fagron	Sodium edetate	13K19-B02-289035	1
Fagron	Sodium edetate	13/23-N03	2
Fagron	Sodium edetate	13I23-N03	1
Fagron	Sodium edetate	16F21-B02-332039	1
Fagron	Sodium edetate	16F21-B02-329900	1
Fagron	Sodium edetate	13K19-B02-290051	2
Fagron	Sodium edetate	16F21-B02	2
Fagron	Sodium edetate	601EC350	1
Fagron	Sodium edetate	16F21B02327696	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Sodium edetate	16f21-b02-327696	1
Fagron	Sodium edetate	13G15-N03	2
Fagron	Sodium edetate	914E-01925	1
Fagron	Sodium edetate	11A31-N06	1
Fagron	Sodium edetate	16F21-B02-323074	1
Fagron	Sodium edetate	16F21-B02-325085	6
Fagron	Sodium edetate	16F21-B02-327696	6
Fagron	Sodium edetate	5012Q-01925	1
Fagron	Sodium edetate	13k19-b02-290051	2
Fagron	Sodium hydrogen carbo...	4413I01892	1
Fagron	Sodium hydrogen carbo...	12G11-N03	7
Fagron	Sodium hydrogen carbo...	12D11-No2	1
Fagron	Sodium hydrogen carbo...	11/12-N16	1
Fagron	Sodium hydrogen carbo...	11I12-N16	3
Fagron	Sodium hydrogen carbo...	12D11-N02	9
Fagron	Sodium hydrogen carbo...	RBA3421100	3
Fagron	Sodium hydrogen carbo...	11K25-N05	3
Fagron	Sodium hydrogen carbo...	13B18-N03	6
Fagron	Sodium hydrogen carbo...	RBA3651200	1
Fagron	Sodium hydrogen carbo...	13D29-N09	4
Fagron	Sodium hydrogen carbo...	13G11-N03	2
Fagron	Sodium hydroxide pell...	12E02-N15	1
Fagron	Sodium hydroxide pell...	13k22-no5	1
Fagron	Sodium hydroxide pell...	13F06-N02	3
Fagron	Sodium monohydrogenph...	13H01-B01	1
Fagron	Sodium monohydrogenph...	13D24-N12	2
Fagron	Sodium monohydrogenph...	613I-01914	1
Fagron	Sodium monohydrogenph...	1513I-01914	1
Fagron	Sodium monohydrogenph...	B44888A	1
Fagron	Sodium monohydrogenph...	2111E-01914	1
Fagron	Sodium monohydrogenph...	11K30-N01	1
Fagron	Sodium salicylate, po...	102091210	1
Fagron	Sodium sulfate decahy...	2134	1
Fagron	Sodium sulfate decahy...	1464	1
Fagron	Sodium tetraborate	12H09-N03	1
Fagron	Sodium thiosulfate pe...	3313U-01922	1
Fagron	Sodium thiosulfate pe...	13092013E	1
Fagron	Sodium thiosulfate pe...	12F20-N10	1
Fagron	Sodium thiosulfate pe...	12F20-N20	1
Fagron	Sorbic acid	121394	1
Fagron	Sorbitan stearate	13D10_N10	1
Fagron	Sorbitan stearate	14J07-B02-299399	1
Fagron	Sorbitan stearate	12E30-N14	1
Fagron	Sorbitan stearate	13D10-N10	1
Fagron	Sorbitan stearate	5052	1
Fagron	Sorbitol	12E30-N24	3
Fagron	Sorbitol	M721140	1
Fagron	Sorbitol	13G04-N04	1
Fagron	Sorbitol	14E21-B02	1
Fagron	Sorbitol	12J15-N01	1
Fagron	Sorbitol	13D08-N27	2
Fagron	Sorbitol	11I09-N11	1
Fagron	Sorbitol	M721140/12L17-N12	1
Fagron	Sorbitol	12L17-N12	3
Fagron	Sorbitol	13E07-B40	2
Fagron	Spironolactone	17E22-B11-335454	2
Fagron	Spironolactone	16I06-B04	1
Fagron	Spironolactone	18A24-B03	1
Fagron	Spironolactone	12J15-N04	1
Fagron	Spironolactone	10H24-N10	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Spironolactone	12/17-N01	1
Fagron	Spironolactone	13G03-N09	1
Fagron	Spironolactone	13G16-B08-288051	1
Fagron	Spironolactone	13D12-N08	1
Fagron	Spironolactone	14B20-B06-294618	1
Fagron	Spironolactone	14B20B06294618	1
Fagron	Spironolactone	14J02-B02-303621	5
Fagron	Spironolactone	14J02-B303621	1
Fagron	Spironolactone	14j02-b02-303621	1
Fagron	Spironolactone	16A14-B05	1
Fagron	Spironolactone	80616-7	1
Fagron	Spironolactone	2415E-03249	1
Fagron	Spironolactone	16A14-B05-317123	1
Fagron	Spironolactone	16I06-B04-329443	1
Fagron	Spironolactone	15G22-B10-311159	2
Fagron	Spironolactone	16D15-B06-320903	2
Fagron	Spironolactone	16A14-B05-525187	1
Fagron	Spironolactone	16E20-B06-325478	2
Fagron	Spironolactone	16E20-B06-322016	2
Fagron	Spironolactone	12/17-n01	1
Fagron	Spironolactone	12I17-N01	1
Fagron	Spironolactone	14B20-B06-294618	4
Fagron	Spironolactone	17H16-B05	2
Fagron	Spironolactone	17H16-B05-344771	3
Fagron	Spironolactone	10C29-N09	1
Fagron	Spironolactone	16I06B04329443	1
Fagron	Sucrose	12G26-N01	2
Fagron	Sucrose	1011208731	2
Fagron	Sucrose	13G04-N01	1
Fagron	Sucrose	2102308731	1
Fagron	Sucrose	13D29-N08	1
Fagron	Sucrose	13G03-N16	1
Fagron	Sucrose	18082014D	1
Fagron	Sucrose	12E29-N04	1
Fagron	Sucrose	12-E03-M01	1
Fagron	Sucrose	11I16-N01	1
Fagron	Sucrose	11/16-N01	1
Fagron	Sucrose	5112I-02207	1
Fagron	Sucrose	11E19-N04	1
Fagron	Sucrose	11/16-N01	1
Fagron	Sucrose	11/20-N05	1
Fagron	Sucrose	11L09-N16	2
Fagron	Sucrose	5071201	1
Fagron	Sucrose	13A07-N03	1
Fagron	Sucrose	13I30-N05	1
Fagron	Sulfanilamide	11L01-No8	1
Fagron	Sweet fennel	2870A130304	3
Fagron	SyrSpend [®] SF Alka	15K03-B01-525173	1
Fagron	SyrSpend [®] SF Alka	18821-b08-347192	1
Fagron	SyrSpend [®] SF Alka	18821-B08-347192	1
Fagron	SyrSpend [®] SF Alka	17D11-804-332593	1
Fagron	SyrSpend [®] SF Alka	18A26-B05-346410	1
Fagron	SyrSpend [®] SF Alka	18A26-Bo5-346410	1
Fagron	SyrSpend [®] SF Alka	17F14-B02-335945	1
Fagron	SyrSpend [®] SF Alka	17D11-bo4-332593	1
Fagron	SyrSpend [®] SF Alka	17D11-B04	2
Fagron	SyrSpend [®] SF Alka	17D11-B04-332593	2
Fagron	SyrSpend [®] SF Alka	18H01-B12-353473	2
Fagron	SyrSpend [®] SF Alka	16B03-B07-525234	1
Fagron	SyrSpend [®] SF pH4	18A22-B05-347899	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	SyrSpend [®] SF pH4	16J10-B03-330757	4
Fagron	SyrSpend [®] SF pH4	18A22-B05-346147	1
Fagron	SyrSpend [®] SF pH4	16J10B03	1
Fagron	SyrSpend [®] SF pH4	16/10-B03	1
Fagron	Tannin	13H30-B03-287645	2
Fagron	Tannin	12H14N05	1
Fagron	Tannin	12H14-N05	6
Fagron	Tannin	12A27-N01	2
Fagron	Tannin	13B28-N08	1
Fagron	Tannin	12E21-N04	4
Fagron	Tannin	11C09-N05	1
Fagron	Tannin	16D08-B04-337161	1
Fagron	Tannin	13H30-B03287645	1
Fagron	Tannin	14E07-B03-298480	2
Fagron	Tannin	14J06-B04-0858	1
Fagron	Tannin	160300126	1
Fagron	Tannin	160300146	2
Fagron	Tannin	2212M-01052	1
Fagron	Tannin	14A07-B01-293468	5
Fagron	Tannin	14J06-B04	1
Fagron	Tannin	14J06-B04-299398	1
Fagron	Tannin	14J06B04299398	1
Fagron	Tannin	14e07-b03-298480	1
Fagron	Tannin	160300029	1
Fagron	Tannin	1014A-01052	1
Fagron	Tannin	15F23-B01-326357	2
Fagron	Tannin	15F23-B01	2
Fagron	Tannin	15F23B01326357	1
Fagron	Tannin	160300178B	2
Fagron	Tannin	16D08-b04	1
Fagron	Tannin	16D08-B04-334834	4
Fagron	Tannin	16D08B04335787	1
Fagron	Tannin	16D08-B04-335787	1
Fagron	Tannin	16D08-B04	1
Fagron	Tannin	16D08B04354470	1
Fagron	Tannin	15F23-B01-322984	2
Fagron	Tartaric acid	12131	2
Fagron	Tartaric acid	12225	1
Fagron	Tartaric acid	13/25-B01-285488	2
Fagron	Tartaric acid	13I25-B01	1
Fagron	Tartaric acid	14G07-B02	1
Fagron	Tartaric acid	1106211053	1
Fagron	Tartaric acid	11E24-N06	1
Fagron	Tartaric acid	14519	1
Fagron	Tartaric acid	14G07-B02-308554	3
Fagron	Tartaric acid	14186212	1
Fagron	Tartaric acid	15045303	2
Fagron	Tartaric acid	1325-B01-294093	1
Fagron	Tartaric acid	13I25-B01-294095	1
Fagron	Tartaric acid	20151112-2	1
Fagron	Tartaric acid	xxxxxxx	1
Fagron	Tartaric acid	1325-B01-294094	1
Fagron	Tartaric acid	12F18-No2	1
Fagron	Tartaric acid	12F18-N02	1
Fagron	Tartaric acid	13I25-B01-294094	3
Fagron	Tartaric acid	2513M01053	1
Fagron	Tartaric acid	12K16-N02	1
Fagron	Tartaric acid	13I25-B01-285487	1
Fagron	Tartaric acid	13B15-N01	10
Fagron	Tartaric acid	12H02-N07	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Tartaric acid	13/25-B01-285487	1
Fagron	Tartaric acid	13B15N01	1
Fagron	Tartaric acid	12K29-B01-285247	1
Fagron	Tartaric acid	11E24-N060752/764	1
Fagron	Tartaric acid	13I25-B01-285488	2
Fagron	Tartaric acid	13I25-B01-285488	1
Fagron	Tartaric acid	13I25-B01-294093	1
Fagron	Tartaric acid	14G07-802	1
Fagron	Tartaric acid	13/25-B01	1
Fagron	Tartaric acid	2414E1053	1
Fagron	Testosterone	13G10-N03	1
Fagron	Testosterone	13g10n03	1
Fagron	Testosterone	1712-B05	1
Fagron	Testosterone	17A12-B05-331496	1
Fagron	Testosterone	13109b03294149	1
Fagron	Testosterone	12F13-N02	1
Fagron	Testosterone	17J26-B01-342033	2
Fagron	Testosterone	13109b03298118	1
Fagron	Testosterone	13G10-N03	1
Fagron	Testosterone	17J26-B01-342637	4
Fagron	Testosterone	13L09-B03-294148	1
Fagron	Testosterone	15E18-B08	1
Fagron	Testosterone	11/07-N12	1
Fagron	Testosterone	GTJ151103	1
Fagron	Testosterone	13I26-N08	1
Fagron	Testosterone propionate		1
Fagron	Testosterone propionate	20121012-3	1
Fagron	Testosterone propionate	12J10-N08	7
Fagron	Testosterone propionate	6081201	1
Fagron	Testosterone propionate	14A28-B12-296960	1
Fagron	Testosterone propionate	12E29-N02	1
Fagron	Testosterone propionate	12E04-N01	2
Fagron	Testosterone propionate	11K18-N04	3
Fagron	Testosterone propionate	11C21-N09	1
Fagron	Testosterone propionate	12A02-N05	3
Fagron	Testosterone propionate	24411M1	1
Fagron	Testosterone propionate	13K25B05287648	1
Fagron	Testosterone propionate	13I03-N02	1
Fagron	Testosterone propionate	13K25-B05-287648	5
Fagron	Testosterone propionate	13F05-N04	14
Fagron	Testosterone propionate	12J10N-08	1
Fagron	Testosterone propionate	32012M1	1
Fagron	Testosterone propionate	3613M-02376	1
Fagron	Testosterone propionate	12/10-N08	1
Fagron	Testosterone propionate	13C08N02	1
Fagron	Testosterone propionate	13K25-B05-291974	2
Fagron	Testosterone propionate	11/12-N03	1
Fagron	Testosterone propionate	36513M1	2
Fagron	Testosterone propionate	12C01-N01	1
Fagron	Testosterone propionate	14A28-B12-291977	1
Fagron	Testosterone propionate	13C08-N02	16
Fagron	Testosterone propionate	13/03-N02	1
Fagron	Testosterone propionate	10L03-N16	3
Fagron	Testosterone propionate	13K25-B05-287647	8
Fagron	Testosterone propionate	15I01-B09-312867	1
Fagron	Testosterone propionate	17B09-B03-332376	2
Fagron	Testosterone propionate	15/01-B09-318371	2
Fagron	Testosterone propionate	15I01-B09-318371	1
Fagron	Testosterone propionate	13K25-B05-291974	1
Fagron	Testosterone propionate	15/01-B09	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Testosterone propionate	18B06-B04-352492	1
Fagron	Testosterone propionate	16E24B01322938	1
Fagron	Testosterone propionate	15I01-B01-318371	1
Fagron	Testosterone propionate	17J16B04346529	1
Fagron	Testosterone propionate	14D15-B03-296959	2
Fagron	Testosterone propionate	16E24-B01-322593	7
Fagron	Testosterone propionate	16C24B04320784	1
Fagron	Testosterone propionate	16E24-B01	1
Fagron	Testosterone propionate	16H22-B05-326483	7
Fagron	Testosterone propionate	16H22-B05-326482	4
Fagron	Testosterone propionate	18/12-B03-355296	1
Fagron	Testosterone propionate	18B06-B04-348557	1
Fagron	Testosterone propionate	18B06-B04-351510	1
Fagron	Testosterone propionate	16E24-B01-325476	1
Fagron	Testosterone propionate	17E22-B09-338592	1
Fagron	Testosterone propionate	16H22B05326483	1
Fagron	Testosterone propionate	17J16-B04-346529	1
Fagron	Testosterone propionate	13K25-B05-287647	1
Fagron	Testosterone propionate	17H25-B02-341321	2
Fagron	Testosterone propionate	16H22-B05	2
Fagron	Testosterone propionate	13F05-N04	1
Fagron	Testosterone propionate	17J16-B03-343010	2
Fagron	Testosterone propionate	13K25-B05-287648	1
Fagron	Testosterone propionate	14L31-B02	1
Fagron	Testosterone propionate	18b14-b02-353456	1
Fagron	Testosterone propionate	18B06-B04-352493	1
Fagron	Testosterone propionate	14A28-B12-296958	2
Fagron	Testosterone propionate	15D03-B02-310694	1
Fagron	Testosterone propionate	14A28-B12	1
Fagron	Testosterone propionate	16C24-B04-320784	5
Fagron	Testosterone propionate	15K03-B03-315218	7
Fagron	Testosterone propionate	18/12-B03-356703	1
Fagron	Testosterone propionate	1501-B09-312867	1
Fagron	Testosterone propionate	14L31B02304722	1
Fagron	Testosterone propionate	14L31-B02-304722	12
Fagron	Testosterone propionate	14L31-B02-304721	8
Fagron	Testosterone propionate	18B06-B04	1
Fagron	Testosterone propionate	14L31B02-304721	1
Fagron	Testosterone propionate	14D15-B03-296959	3
Fagron	Testosterone propionate	14A28-B12	1
Fagron	Testosterone propionate	15K03-B-03-315218	1
Fagron	Testosterone propionate	11G15-N07	1
Fagron	Testosterone propionate	36814M1	1
Fagron	Testosterone propionate	15k03-b03	1
Fagron	Testosterone propionate	14L31-B02	1
Fagron	Tetracaine base	14a16-b04-290208	1
Fagron	Tetracaine base	15/16-B06-316260	1
Fagron	Tetracaine base	15/16-B06-312783	2
Fagron	Tetracaine base	14G30-B08	1
Fagron	Tetracaine base	15G30-B04-312060	1
Fagron	Tetracaine base	14g30-b08-297306	1
Fagron	Tetracaine base	1516-B06-316259	1
Fagron	Tetracaine base	12F29-N05	7
Fagron	Tetracaine base	1516-B06-316260	1
Fagron	Tetracaine base	14A16-B04	1
Fagron	Tetracaine base	2913Q-07165	1
Fagron	Tetracaine base	11D07-N20	1
Fagron	Tetracaine base	14G30-B08-297307	1
Fagron	Tetracaine base	16E23-B04-323269	1
Fagron	Tetracaine base	14A16-B04-290209	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Tetracaine base	13D30-N11	2
Fagron	Tetracaine base	11101-N09	1
Fagron	Tetracaine base	12F26-N05	1
Fagron	Tetracaine base	14G30-B08	1
Fagron	Tetracaine base	14K03-B02	1
Fagron	Tetracaine base	14g30-b08-297305	1
Fagron	Tetracaine base	11L09-N03	3
Fagron	Tetracaine base	12L19-N10	5
Fagron	Tetracaine base	14G30-B08-297306	3
Fagron	Tetracaine base	TC0090913	1
Fagron	Tetracaine base	14A16-B04-2490207	1
Fagron	Tetracaine base	11G12-N13	3
Fagron	Tetracaine base	15.16-B06-312783	1
Fagron	Tetracaine base	12H23-N07	1
Fagron	Tetracaine base	12f29-n05	1
Fagron	Tetracaine base	5111202	1
Fagron	Tetracaine base	TC0090913	1
Fagron	Tetracaine base	14K03B02303178	1
Fagron	Tetracaine base	14G30-B08-297305	1
Fagron	Tetracaine base	16J12-B06	1
Fagron	Tetracaine base	16G29B02325570	2
Fagron	Tetracaine base	17K13-B02-344448	3
Fagron	Tetracaine base	17K13-B02-347172	1
Fagron	Tetracaine base	17E18-B02-337216	1
Fagron	Tetracaine base	17E18-B02-337217	2
Fagron	Tetracaine base	14K03-B02-303179	2
Fagron	Tetracaine base	16J12-B06-329087	1
Fagron	Tetracaine base	16A13-B07-322947	1
Fagron	Tetracaine base	16G29-B02-325570	1
Fagron	Tetracaine base	16J12-B06-327371	2
Fagron	Tetracaine hydrochlor...	12K01-N04	7
Fagron	Tetracaine hydrochlor...	10A29-N12	2
Fagron	Tetracaine hydrochlor...	12H24-N05	15
Fagron	Tetracaine hydrochlor...	11C22-N04	1
Fagron	Tetracaine hydrochlor...	11A31-N15	1
Fagron	Tetracaine hydrochlor...	12C22-N04	5
Fagron	Tetracaine hydrochlor...	3913I-02039	1
Fagron	Tetracaine hydrochlor...	3413U02039	1
Fagron	Tetracaine hydrochlor...	13D11-N07	12
Fagron	Tetracaine hydrochlor...	12H-N05	1
Fagron	Tetracaine hydrochlor...	11/15-N12	1
Fagron	Tetracaine hydrochlor...	12H24-N050696/1050	1
Fagron	Tetracaine hydrochlor...	13D11-N01	1
Fagron	Tetracaine hydrochlor...	212131/221630-0002	1
Fagron	Tetracaine hydrochlor...	17G03-B07-351509	1
Fagron	Tetracaine hydrochlor...	17G03-B07-349368	1
Fagron	Tetracaine hydrochlor...	17G03-B07-345008	4
Fagron	Tetracaine hydrochlor...	17Go3-B07-341853	1
Fagron	Tetracaine hydrochlor...	13G02-N01	1
Fagron	Tetracaine hydrochlor...	17114181	2
Fagron	Tetracaine hydrochlor...	13H21-B02-286558	4
Fagron	Tetracaine hydrochlor...	13H21-B02	2
Fagron	Tetracaine hydrochlor...	17G03-B07-338540	1
Fagron	Tetracaine hydrochlor...	17G03-B07-338539	2
Fagron	Tetracaine hydrochlor...	4313E-02039	1
Fagron	Tetracaine hydrochlor...	13K25-B04-290278	1
Fagron	Tetracaine hydrochlor...	17G03-B07-341319	1
Fagron	Tetracaine hydrochlor...	10990841	2
Fagron	Tetracaine hydrochlor...	15K04B03318369	1
Fagron	Tetracaine hydrochlor...	15K04-B03-318368	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Tetracaine hydrochlor...	15H24-B03	3
Fagron	Tetracaine hydrochlor...	13K25B04290277	1
Fagron	Tetracaine hydrochlor...	14C18-B03	1
Fagron	Tetracaine hydrochlor...	15K04-bo3-312151	1
Fagron	Tetracaine hydrochlor...	13K25-B04	1
Fagron	Tetracaine hydrochlor...	16G28-B03-327369	2
Fagron	Tetracaine hydrochlor...	17G05-B04-341204	1
Fagron	Tetracaine hydrochlor...	10990834	3
Fagron	Tetracaine hydrochlor...	14C18-Bo3-300200	1
Fagron	Tetracaine hydrochlor...	14C18-B03-300199	5
Fagron	Tetracaine hydrochlor...	16G28-B03-327368	15
Fagron	Tetracaine hydrochlor...	17G03-B07	5
Fagron	Tetracaine hydrochlor...	16G2-B03-327368	1
Fagron	Tetracaine hydrochlor...	218I-02039	1
Fagron	Tetracaine hydrochlor...	17G03-B07-338538	6
Fagron	Tetracaine hydrochlor...	16G28-B03-334124	2
Fagron	Tetracaine hydrochlor...	14710369	2
Fagron	Tetracaine hydrochlor...	16G28-B03-330363	1
Fagron	Tetracaine hydrochlor...	1418-03-300199	1
Fagron	Tetracaine hydrochlor...	15K04-B03_322151	1
Fagron	Tetracaine hydrochlor...	15H24-B03-313594	5
Fagron	Tetracaine hydrochlor...	14C18-B03-293966	7
Fagron	Tetracaine hydrochlor...	14c18d03300199	1
Fagron	Tetracaine hydrochlor...	10990837	3
Fagron	Tetracaine hydrochlor...	14G14-B08-300198	6
Fagron	Tetracaine hydrochlor...	13K25-B04290276	1
Fagron	Tetracaine hydrochlor...	13K25-B04-290276	3
Fagron	Tetracaine hydrochlor...	10330837	1
Fagron	Tetracaine hydrochlor...	15A12-B03-308857	8
Fagron	Tetracaine hydrochlor...	16G8-B03-327368	1
Fagron	Tetracaine hydrochlor...	13K25-B04...	1
Fagron	Tetracaine hydrochlor...	15h24-b03-313594	1
Fagron	Tetracaine hydrochlor...	13k25-b04-290276	1
Fagron	Tetracaine hydrochlor...	13K25-Bo4-290276	1
Fagron	Tetracycline hydrochl...	13L10-B02-294454	1
Fagron	Tetracycline hydrochl...	14H06-B05	1
Fagron	Tetracycline hydrochl...	14H06-B05-303911	1
Fagron	Tetracycline hydrochl...	12E10-N05	1
Fagron	Tetracycline hydrochl...	14H06-B05-303912	6
Fagron	Tetracycline hydrochl...	13L10-B02-294455	2
Fagron	Tetracycline hydrochl...	12J15-N05	1
Fagron	Tetracycline hydrochl...	13F17-N03	5
Fagron	Tetracycline hydrochl...	14H06-B05-299498	1
Fagron	Tetracycline hydrochl...	14H06-Bo5-303912	1
Fagron	Tetracycline hydrochl...	13J22-N04	1
Fagron	Tetracycline hydrochl...	14H06-B05-314870	1
Fagron	Tetracycline hydrochl...	14H06-B05-311707	1
Fagron	Tetracycline hydrochl...	18308305	1
Fagron	Thalidomide	13K20-N01	1
Fagron	Thalidomide	14E27-B04-299297	1
Fagron	Thiamine chloride hyd...	13E22-No8	1
Fagron	Thiamine chloride hyd...	10111401	1
Fagron	Thiamine chloride hyd...	17F23-B02-340810	1
Fagron	Thiamine chloride hyd...	190115006	1
Fagron	Thiomersal	12117-N02	1
Fagron	Thiomersal	12J17-N02	1
Fagron	Thiomersal	12F19-N11	1
Fagron	Thymol	12K01-N02	2
Fagron	Thymol	12H29-N02	2
Fagron	Thymol	713575	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Thymol	12M29-N02	1
Fagron	Thymol	12E30-N11	1
Fagron	Thymol	12B06-N11	1
Fagron	Tobramycin sulfate	16B05-B03-317988	2
Fagron	Tobramycin sulfate	17E05-B03-338122	2
Fagron	Tobramycin sulfate	16B05-B03	4
Fagron	Tobramycin sulfate	16B05-B03-330974	2
Fagron	Tobramycin sulfate	13H09-N09	1
Fagron	Tobramycin sulfate	16B05-B03-319433	1
Fagron	Tobramycin sulfate	15/17-B04-315697	1
Fagron	Tosylchloramide sodium	14L12-B01-301915	1
Fagron	Tosylchloramide sodium	11L08-N04	1
Fagron	Tosylchloramide sodium	14B27-B03-291967	1
Fagron	Tragacanth, powdered	13204705	5
Fagron	Tranexamic acid	13D08-N24	1
Fagron	Tranexamic acid	18F20-B01-351693	1
Fagron	Tranexamic acid	5041804	1
Fagron	Tranexamic acid	18F20-B01-351691	2
Fagron	Tranexamic acid	18A17-B05-345038	1
Fagron	Tranexamic acid	13D08-N240002/40	1
Fagron	Tranexamic acid	16B03-B02-319617	2
Fagron	Tranexamic acid	16J29-B04-312372	1
Fagron	Tranexamic acid	15J29-B04-312372	1
Fagron	Tranexamic acid	12D11-N02	1
Fagron	Tranexamic acid	13G24-N21	1
Fagron	Tranexamic acid	13D08-N24	1
Fagron	Tranexamic acid	13K08-B04-287966	1
Fagron	Tranexamic acid	15J29-B04	2
Fagron	Tretinoin	1410001594	1
Fagron	Tretinoin	20120601	2
Fagron	Tretinoin	13D30-N14-3365319091	1
Fagron	Tretinoin	14L24-B01-303152	1
Fagron	Tretinoin	14D09-B04-295551	1
Fagron	Tretinoin	13B05-B02-287360	5
Fagron	Tretinoin	13B05-B02-287361	2
Fagron	Tretinoin	1440104	1
Fagron	Tretinoin	13C11-N06	2
Fagron	Tretinoin	14D09-B04-293594	1
Fagron	Tretinoin	13D30N14	1
Fagron	Tretinoin	14093505	1
Fagron	Tretinoin	11K29-N01	1
Fagron	Tretinoin	12106-N01	1
Fagron	Tretinoin	14L24-B01	1
Fagron	Tretinoin	14D09-B04-293225	2
Fagron	Tretinoin	231115T	1
Fagron	Tretinoin	15F01-B06-308856	3
Fagron	Tretinoin	14L24-B01-303157	5
Fagron	Tretinoin	16F14B04325560	1
Fagron	Tretinoin	16F14-B04-326181	3
Fagron	Tretinoin	16K03-B06-328388	5
Fagron	Tretinoin	17K01B05349946	1
Fagron	Tretinoin	17K01-B05-347436	1
Fagron	Tretinoin	17K01-B05	2
Fagron	Tretinoin	17K01-B05-345784	2
Fagron	Tretinoin	17357007	6
Fagron	Tretinoin	18A10-B03-344364	1
Fagron	Tretinoin	16F14-B04-325567	1
Fagron	Tretinoin	17K01-B05-341941	1
Fagron	Tretinoin	16K03-B06	3
Fagron	Tretinoin	170706TTRER-20151004	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Tretinoin	17K01-B05-341936	3
Fagron	Tretinoin	16F14-B04-325560	2
Fagron	Tretinoin	16F14-B04	4
Fagron	Tretinoin	17K01-B05-347435	1
Fagron	Tretinoin	15F01-B06-308722	6
Fagron	Tretinoin	15F01-B06-317884	2
Fagron	Tretinoin	14009-804	1
Fagron	Triamcinolone acetonide	13121201	1
Fagron	Triamcinolone acetonide	13k14-B04-289032	1
Fagron	Triamcinolone acetonide	9041305	1
Fagron	Triamcinolone acetonide	13K14-B04-289033	10
Fagron	Triamcinolone acetonide	13K14-B04-289033	1
Fagron	Triamcinolone acetonide	13k15-b31-289030	1
Fagron	Triamcinolone acetonide	13i26-N11	1
Fagron	Triamcinolone acetonide	14F30-B01-296287	3
Fagron	Triamcinolone acetonide	NT1430M	4
Fagron	Triamcinolone acetonide	14G23-B07-297620	2
Fagron	Triamcinolone acetonide	13K14-B04	2
Fagron	Triamcinolone acetonide	3212U-02516 (NP)	1
Fagron	Triamcinolone acetonide	13C13-20	1
Fagron	Triamcinolone acetonide	NS2236M	1
Fagron	Triamcinolone acetonide	NT1094M	1
Fagron	Triamcinolone acetonide	5061303	1
Fagron	Triamcinolone acetonide	30.06.2017	1
Fagron	Triamcinolone acetonide	13K14-B04-289032	4
Fagron	Triamcinolone acetonide	1210-N09	1
Fagron	Triamcinolone acetonide	13K15-B31-289031	5
Fagron	Triamcinolone acetonide	1326-N11	1
Fagron	Triamcinolone acetonide	NS0571M	1
Fagron	Triamcinolone acetonide	13126-N11	1
Fagron	Triamcinolone acetonide	11B22-N05	1
Fagron	Triamcinolone acetonide	12E15-N03	23
Fagron	Triamcinolone acetonide	C22H25CIN209	1
Fagron	Triamcinolone acetonide	12I10-N09	6
Fagron	Triamcinolone acetonide	12I10-N09	3
Fagron	Triamcinolone acetonide	NU0930M	1
Fagron	Triamcinolone acetonide	5061205	1
Fagron	Triamcinolone acetonide	14G23-B07-297619	2
Fagron	Triamcinolone acetonide	12B08-N01	7
Fagron	Triamcinolone acetonide	12/10-N09	7
Fagron	Triamcinolone acetonide	14G23-B07-297620	1
Fagron	Triamcinolone acetonide	11K11-N10	4
Fagron	Triamcinolone acetonide	120600801	1
Fagron	Triamcinolone acetonide	14F30-B01	1
Fagron	Triamcinolone acetonide	NS0996M	2
Fagron	Triamcinolone acetonide	14G23-B07-297620	1
Fagron	Triamcinolone acetonide	14G23-B07-297618	3
Fagron	Triamcinolone acetonide	20031302	1
Fagron	Triamcinolone acetonide	11N11-N10	1
Fagron	Triamcinolone acetonide	14K24N-N01	1
Fagron	Triamcinolone acetonide	18C15-B05-350263	1
Fagron	Triamcinolone acetonide	13K15-B31	5
Fagron	Triamcinolone acetonide	16D14-B11-328167	11
Fagron	Triamcinolone acetonide	16D04-B11-325479	1
Fagron	Triamcinolone acetonide	16D14_B11_328167	1
Fagron	Triamcinolone acetonide	16K07-B06-329445	3
Fagron	Triamcinolone acetonide	NW2261M	1
Fagron	Triamcinolone acetonide	16B04-B07318664	1
Fagron	Triamcinolone acetonide	18c15-b05-349087	1
Fagron	Triamcinolone acetonide	16K07-B06-329446	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Triamcinolone acetonide	16D14-B11-325480	3
Fagron	Triamcinolone acetonide	NT0770M	1
Fagron	Triamcinolone acetonide	16K07-B06-332101	6
Fagron	Triamcinolone acetonide	16B04-Bo7-318664	1
Fagron	Triamcinolone acetonide	NT0487M	1
Fagron	Triamcinolone acetonide	16D14-B11-325481	1
Fagron	Triamcinolone acetonide	16B04	1
Fagron	Triamcinolone acetonide	16014-B12-3231318	1
Fagron	Triamcinolone acetonide	13K15-B31-289025	1
Fagron	Triamcinolone acetonide	16D14-B11	3
Fagron	Triamcinolone acetonide	16D14-B11-325479	6
Fagron	Triamcinolone acetonide	1810015-01	14
Fagron	Triamcinolone acetonide	16D14-B11-322135	4
Fagron	Triamcinolone acetonide	15G31B07313597	1
Fagron	Triamcinolone acetonide	13I26-N11	12
Fagron	Triamcinolone acetonide	16D14-B11-325482	4
Fagron	Triamcinolone acetonide	16D14-B11_322158	1
Fagron	Triamcinolone acetonide	16D14B11325479	1
Fagron	Triamcinolone acetonide	16-B07-318664	1
Fagron	Triamcinolone acetonide	17H03-B12-343457	4
Fagron	Triamcinolone acetonide	17c09-b01-340174	1
Fagron	Triamcinolone acetonide	17H03-B12-343458	3
Fagron	Triamcinolone acetonide	17C-09-B01-334197	1
Fagron	Triamcinolone acetonide	17H03-B12-343672	5
Fagron	Triamcinolone acetonide	13K15B31289029	1
Fagron	Triamcinolone acetonide	NT077oM	1
Fagron	Triamcinolone acetonide	17H03-B12-343459	1
Fagron	Triamcinolone acetonide	17C09-B01-340173	3
Fagron	Triamcinolone acetonide	17C09-B01-334197	3
Fagron	Triamcinolone acetonide	NX1540M	1
Fagron	Triamcinolone acetonide	18C15-B05	3
Fagron	Triamcinolone acetonide	18C15-B05-347275	4
Fagron	Triamcinolone acetonide	17H03-B12	4
Fagron	Triamcinolone acetonide	18C15-B05_348063	1
Fagron	Triamcinolone acetonide	17H03-B12-345552	3
Fagron	Triamcinolone acetonide	18C15-B05-348063	2
Fagron	Triamcinolone acetonide	18C15B05347275	2
Fagron	Triamcinolone acetonide	17C09-B01	2
Fagron	Triamcinolone acetonide	18C15-B05-349087	3
Fagron	Triamcinolone acetonide	ny0485m	1
Fagron	Triamcinolone acetonide	17C09-B01_337464	1
Fagron	Triamcinolone acetonide	NX0241M	1
Fagron	Triamcinolone acetonide	17C09-B01-334196	3
Fagron	Triamcinolone acetonide	17C09-B01-337464	2
Fagron	Triamcinolone acetonide	16D14-B11-322158	7
Fagron	Triamcinolone acetonide	16D14-B11-321318	3
Fagron	Triamcinolone acetonide	17C09-B01-340174	3
Fagron	Triamcinolone acetonide	18F14-B01-352343	1
Fagron	Triamcinolone acetonide	15G31-B07	1
Fagron	Triamcinolone acetonide	15G31-B07-314045	6
Fagron	Triamcinolone acetonide	15G31-B07-313597	5
Fagron	Triamcinolone acetonide	1531-07-313596	1
Fagron	Triamcinolone acetonide	15G31-B07-313596	6
Fagron	Triamcinolone acetonide	15G31-B07-314047	2
Fagron	Triamcinolone acetonide	15C19-B06-309194	2
Fagron	Triamcinolone acetonide	16B04-B07-319471	4
Fagron	Triamcinolone acetonide	15B02-B03-303828	1
Fagron	Triamcinolone acetonide	14F30-B021-296287	1
Fagron	Triamcinolone acetonide	14355709	4
Fagron	Triamcinolone acetonide	NV0262M	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Triamcinolone acetonide	15C19-B06-308571	1
Fagron	Triamcinolone acetonide	15C19-B06-308572	7
Fagron	Triamcinolone acetonide	Nw0635M	1
Fagron	Triamcinolone acetonide	15B02-B03-303830	1
Fagron	Triamcinolone acetonide	14K14-B05-302887	2
Fagron	Triamcinolone acetonide	15B02-B03	4
Fagron	Triamcinolone acetonide	17C09-B01-336881	1
Fagron	Triamcinolone acetonide	15g31-b07-314046	1
Fagron	Triamcinolone acetonide	15G31B07	2
Fagron	Triamcinolone acetonide	13c13-n20	1
Fagron	Triamcinolone acetonide	15B02-B03-303828	1
Fagron	Triamcinolone acetonide	11K11-N10	1
Fagron	Triamcinolone acetonide	13e15-n03	1
Fagron	Triamcinolone acetonide	17h03-b12-343672	1
Fagron	Triamcinolone acetonide	13C13-N20	3
Fagron	Triamcinolone acetonide	411082-1-P	1
Fagron	Triamcinolone acetonide	13E15-N03	24
Fagron	Triamcinolone acetonide	16B04-B07-318664	4
Fagron	Triamcinolone acetonide	18K16-B05	1
Fagron	Triamcinolone acetonide	15G31-B07-314046	2
Fagron	Triamcinolone acetonide		1
Fagron	Triamcinolone acetonide	12E15-N03	1
Fagron	Triamcinolone acetonide	16B04-B07	4
Fagron	Triamcinolone acetonide	12B08N01	1
Fagron	Triamcinolone acetonide	12/-N09	1
Fagron	Triamcinolone acetonide	15L21-B03	1
Fagron	Triamcinolone acetonide	15C19-B06	1
Fagron	Triamcinolone acetonide	16B04-B07-318957	6
Fagron	Triamcinolone acetonide	13E15N03	1
Fagron	Triamcinolone acetonide	13E15N03	2
Fagron	Triamcinolone acetoni...	12I07-No4	1
Fagron	Triamcinolone acetoni...	13A15-N05	1
Fagron	Triamcinolone acetoni...	12A05-N07	3
Fagron	Triamcinolone acetoni...	13F21-N01	2
Fagron	Triamcinolone acetoni...	11G14-N02	1
Fagron	Triamcinolone acetoni...	12I07-N04	1
Fagron	Triamcinolone acetoni...	12I07-N04	1
Fagron	Triamcinolone acetoni...	12G23-N06	2
Fagron	Triamcinolone acetoni...	13F20-N010113/150	1
Fagron	Triamcinolone acetoni...	12I07-N04022/359	1
Fagron	Triamcinolone acetoni...	12F21-N01	1
Fagron	Triamcinolone acetoni...	12E23-N01	2
Fagron	Triamcinolone acetoni...	14F17-N02	1
Fagron	Triamcinolone, micron...	13C26-N04	1
Fagron	Triamcinolone, micron...	12Do2-N14	1
Fagron	Triamcinolone, micron...	12F01-N19	2
Fagron	Triamcinolone, micron...	13D30-N15	2
Fagron	Triamcinolone, micron...	12D02-N14	1
Fagron	Triamcinolone, micron...	14B13-B30-316102	1
Fagron	Triamcinolone, micron...	14B13-B30-312868	1
Fagron	Triamcinolone, micron...	15B13-B30-312868	1
Fagron	Triamcinolone, micron...	12A03-N09	1
Fagron	Triamcinolone, micron...	NT1484M	1
Fagron	Triamcinolone, micron...	14B13-B30	1
Fagron	Triamcinolone, micron...	11G26-N04	1
Fagron	Triamcinolone, micron...	11D28-N04	1
Fagron	Triamcinolone, micron...	16A13-B06	1
Fagron	Triamcinolone, micron...	14B13-B30-293473	1
Fagron	Triamcinolone, micron...	18A19-B05	1
Fagron	Trichloroacetic acid	13B19-M04	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Trichloroacetic acid	13H22-M04	1
Fagron	Triclosan	15G07-B03-310041	1
Fagron	Triclosan	14F23-B02_295756	1
Fagron	Triclosan	13L18-B35_289995	1
Fagron	Triclosan	13L-18-B35-289985	1
Fagron	Triclosan	13L18-B35-289994	6
Fagron	Triclosan	14F23-B02-300219	9
Fagron	Triclosan	12C21-N05	22
Fagron	Triclosan	13L18-B35-289985	1
Fagron	Triclosan	13L18B35289984	1
Fagron	Triclosan	13H09N01	2
Fagron	Triclosan	13L18-B35-289993	3
Fagron	Triclosan	18082014B	1
Fagron	Triclosan	700228-0005	1
Fagron	Triclosan	Int.Ch.1050814	1
Fagron	Triclosan	12G09-N09	8
Fagron	Triclosan	2914E-03402	1
Fagron	Triclosan	13L18-B35-289987	8
Fagron	Triclosan	13L18-B35-289985	14
Fagron	Triclosan	13H09-No1	1
Fagron	Triclosan	106666066	1
Fagron	Triclosan	10666066	2
Fagron	Triclosan	10112014C	1
Fagron	Triclosan	14F23-B02-295762	3
Fagron	Triclosan	14F23-B02-295760	11
Fagron	Triclosan		2
Fagron	Triclosan	14F23-B02-295756	2
Fagron	Triclosan	10J18-N01	2
Fagron	Triclosan	11H09-N02	6
Fagron	Triclosan	200315B	1
Fagron	Triclosan	4394083	1
Fagron	Triclosan	11238004	2
Fagron	Triclosan	C12H7C1302	1
Fagron	Triclosan	14F23-B02-295761	11
Fagron	Triclosan	11L19-N11	6
Fagron	Triclosan	11J28-N02	5
Fagron	Triclosan	14F-B02-293859	1
Fagron	Triclosan	13L18-35-289995	1
Fagron	Triclosan	14F23-B02-295765	7
Fagron	Triclosan	14F23-B02-293852	2
Fagron	Triclosan	14F23-B02	6
Fagron	Triclosan	14F23-B02-295759	9
Fagron	Triclosan	160614	1
Fagron	Triclosan	13L18-B35-289983	10
Fagron	Triclosan	13H08-N100450/1231	1
Fagron	Triclosan	13092012C	1
Fagron	Triclosan	13Do2-N02	1
Fagron	Triclosan	12G17-N04	38
Fagron	Triclosan	434083	1
Fagron	Triclosan	12G17-n04	1
Fagron	Triclosan	13D02N02	1
Fagron	Triclosan	12G09-no9	1
Fagron	Triclosan	2901313402	1
Fagron	Triclosan	2513Q-03402	1
Fagron	Triclosan	13A11N03	1
Fagron	Triclosan	12G17.N04	1
Fagron	Triclosan	13A11-N03	17
Fagron	Triclosan	13a11n03	1
Fagron	Triclosan	11A05-N13	1
Fagron	Triclosan	12G17N04	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Triclosan	11D20-N10	1
Fagron	Triclosan	12g17-n04	2
Fagron	Triclosan	13H08-N10	34
Fagron	Triclosan	90000383	1
Fagron	Triclosan	13d02-N02	1
Fagron	Triclosan	13H08N10	2
Fagron	Triclosan	13L18-B35-289	1
Fagron	Triclosan	13L18B35289987	1
Fagron	Triclosan	13L18-B35	14
Fagron	Triclosan	13L18-B35-289986	7
Fagron	Triclosan	13L18-B35-289984	7
Fagron	Triclosan	13L18-B35-289988	1
Fagron	Triclosan	13H09-N011296/2856	1
Fagron	Triclosan	1214I03402	1
Fagron	Triclosan	14F23-B02-295757	4
Fagron	Triclosan	10006226	3
Fagron	Triclosan	14022014C	1
Fagron	Triclosan	9000383	2
Fagron	Triclosan	13D02-N-02	1
Fagron	Triclosan	13H08-N-10	1
Fagron	Triclosan	1209213402	1
Fagron	Triclosan	6012014E	1
Fagron	Triclosan	1006226	1
Fagron	Triclosan	13H09-N01	69
Fagron	Triclosan	13h08-n10	1
Fagron	Triclosan	14F23-B02-295763	3
Fagron	Triclosan	12L06-N01	44
Fagron	Triclosan	17G07-B03-310038	1
Fagron	Triclosan	15C10-806-312981	1
Fagron	Triclosan	151010-B06-312984	1
Fagron	Triclosan	190615A	1
Fagron	Triclosan	15g07-b03-310035	1
Fagron	Triclosan	15C10-B06-312980	5
Fagron	Triclosan	15C10-B06-311161	12
Fagron	Triclosan	14F23-B02-295758	10
Fagron	Triclosan	15C10-B06-310032	5
Fagron	Triclosan	14F23B02295758	1
Fagron	Triclosan	15C10B06312983	1
Fagron	Triclosan	15C10-B06-310031	8
Fagron	Triclosan	15C10B06312982	1
Fagron	Triclosan	14F23-B02-2959947	1
Fagron	Triclosan	15C10-B06-312981	25
Fagron	Triclosan	14F23-B02295757	1
Fagron	Triclosan	14F23B02295763	1
Fagron	Triclosan	15G07B03310040	1
Fagron	Triclosan	15C10B06311161	2
Fagron	Triclosan	15C10-B06-312983	7
Fagron	Triclosan	15c10-B06-311161	1
Fagron	Triclosan	15C10-B06-312982	19
Fagron	Triclosan	15G07	1
Fagron	Triclosan	15C10-B06	7
Fagron	Triclosan	15C10-B06_312980	1
Fagron	Triclosan	15c10-b06-312980	1
Fagron	Triclosan	14F23-Bo2-295947	1
Fagron	Triclosan	15C10-B06-3161	1
Fagron	Triclosan	15C-B06-312982	1
Fagron	Triclosan	15G07-B03-310036	4
Fagron	Triclosan	17K02-B06-343544	2
Fagron	Triclosan	15G07-B03	3
Fagron	Triclosan	15G07B03310037	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Triclosan	14f23-b02-295757	1
Fagron	Triclosan	15G07-B03-310038	7
Fagron	Triclosan	17K02-B06-343543	1
Fagron	Triclosan	15G01-B03-310038	1
Fagron	Triclosan	TCS0214037	3
Fagron	Triclosan	15607-B03-310035	1
Fagron	Triclosan	15G07-B03-310039	2
Fagron	Triclosan	17K02-B06-343548	3
Fagron	Triclosan	14F23B02300219	3
Fagron	Triclosan	14F23-B02-295947	5
Fagron	Triclosan	17148812	1
Fagron	Triclosan	17K02-B06	6
Fagron	Triclosan	17081702	2
Fagron	Triclosan	15C10B06310031	2
Fagron	Triclosan	15c10-b06-312981	1
Fagron	Triclosan	K001.1/14/12/823	1
Fagron	Triclosan	K001-1/14/12/823	1
Fagron	Triclosan	12.11.201550g24,15	1
Fagron	Triclosan	14f23-b02-295758	1
Fagron	Triclosan	15G047-B03-0001	1
Fagron	Triclosan	15G07-B03-310037	15
Fagron	Triclosan	15G07-B03-310040	8
Fagron	Trometamol	12I27-No3	1
Fagron	Trometamol	12k14no5	1
Fagron	Trometamol	294101	1
Fagron	Trometamol	13C08-N07	3
Fagron	Trometamol	12K14-N05	2
Fagron	Trometamol	610I-02519	1
Fagron	Trometamol	13E22-N05	6
Fagron	Trometamol	18A23-B10-349673	1
Fagron	Trometamol	12E25-B40	1
Fagron	Trometamol	12I27-N03	1
Fagron	Trometamol	17A16-B04	1
Fagron	Trometamol	17A16-B04-334121	2
Fagron	Trometamol	17A16-B04-340796	1
Fagron	Trometamol	16D14-B09	1
Fagron	Trometamol	14F19-B03	1
Fagron	Trometamol	16D14-B09-328729	2
Fagron	Trometamol	13119-B02	1
Fagron	Trometamol	13/19-B02-293595	2
Fagron	Trometamol	13I19-B02-286158	1
Fagron	Trometamol	13/19-B02	1
Fagron	Trometamol	8386E050	1
Fagron	Trometamol	13119-B02-286158	1
Fagron	Trometamol	13/19-B02-286158	4
Fagron	Trometamol	13I19-B02	1
Fagron	Trometamol	14F19-B03-304293	3
Fagron	Trometamol	16(18A23-B10-347884	1
Fagron	Tropicamide	15F18-B04-308977	1
Fagron	Tropicamide	14104-b02-315694	1
Fagron	Tropicamide	23071508	1
Fagron	Tropicamide	14D23B04293294	1
Fagron	Tropicamide	12104-n03	1
Fagron	Tropicamide	12F18-N06	4
Fagron	Tropicamide	20021302	1
Fagron	Tyrothricin	15A06-B02-303893	1
Fagron	Tyrothricin	17E09-B03-341231	1
Fagron	Tyrothricin	15A06-B02-323005	5
Fagron	Tyrothricin	15a06-B02-323005	1
Fagron	Tyrothricin	13A14-N04	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Tyrothricin	17E09-B03	1
Fagron	Tyrothricin	15A06-B02-3005	1
Fagron	Tyrothricin	17E09-B03-352743	2
Fagron	Tyrothricin	17E09-B03-351982	1
Fagron	Ubidecarenone	14H26-B04-298244	1
Fagron	Ubidecarenone	16B10-B05-327606	1
Fagron	Ubidecarenone	16G26-B01	1
Fagron	Ubidecarenone	14E12-B06-293861	1
Fagron	Ubidecarenone	18a16-b03-351350	1
Fagron	Ubidecarenone	18A16-B03-351350	1
Fagron	Ubidecarenone	16G26-B01-344334	2
Fagron	Urea	14/18-B03-303576	2
Fagron	Urea	14E07B05293469	1
Fagron	Urea	14E07-B05-297165	3
Fagron	Urea	2014162937	3
Fagron	Urea	1418B03303581	1
Fagron	Urea	14I18-B03	5
Fagron	Urea	12F25-N05	5
Fagron	Urea	13/17-1302-287602	1
Fagron	Urea	1211822	1
Fagron	Urea	14118-B03-0303581	1
Fagron	Urea	14I18-B03-303576	5
Fagron	Urea	14E07-B05-297167	6
Fagron	Urea	1418b03303581	1
Fagron	Urea	14I18-B03-303581	4
Fagron	Urea	14E07-B07-293469	1
Fagron	Urea	13D03-N01	9
Fagron	Urea	3613M-02544	1
Fagron	Urea	235/12I04-N06	1
Fagron	Urea	12I04-N06	1
Fagron	Urea	13E17-N02	1
Fagron	Urea	13J17-B02-287604	7
Fagron	Urea	12110-n08	1
Fagron	Urea	13J17-B02-287602	4
Fagron	Urea	13J17-B02	2
Fagron	Urea	12I04-N06	1
Fagron	Urea	13L20-B03-289708	5
Fagron	Urea	14118-B03	1
Fagron	Urea	13L20-B03_289708	1
Fagron	Urea	11H01-N08	2
Fagron	Urea	11C25-N01	1
Fagron	Urea	14E07-B05-293469	13
Fagron	Urea	13L20B03	1
Fagron	Urea	14E07-B05	8
Fagron	Urea	14E-07-B05-293469	1
Fagron	Urea	12C20-N05	1
Fagron	Urea	14E07B05	3
Fagron	Urea	14E8900-0001	1
Fagron	Urea	11F24-N16	1
Fagron	Urea	13J17B02-287604	1
Fagron	Urea	13D03-N001	1
Fagron	Urea	16I13-B02	1
Fagron	Urea	14/18B03303581	1
Fagron	Urea	7L17-N01	1
Fagron	Urea	16/13-B02	1
Fagron	Urea	16113-B02-332953	1
Fagron	Urea	16/13-B02-332953	1
Fagron	Urea	17H16-B01-342080	2
Fagron	Urea	16I13-B02-328152	2
Fagron	Urea	1613b02332953	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Urea	17H16-B01	1
Fagron	Urea	16I13-B02-330166	2
Fagron	Urea	18A24-B07-349015	2
Fagron	Urea	2018163455	1
Fagron	Urea	18A24-B07-352707	1
Fagron	Urea	18A24-B07-347410	1
Fagron	Urea	18G04-B03-35339	1
Fagron	Urea	18G04-B03-353339	3
Fagron	Urea	18A24-B07-351975	3
Fagron	Urea	18G04-B03-355378	1
Fagron	Urea	18G04-B03	1
Fagron	Urea	18A24-B07-346692	1
Fagron	Urea	16D22-B02-326481	8
Fagron	Urea	16/13-B02-330166	1
Fagron	Urea	16D22-802-326481	1
Fagron	Urea	B03-303581	1
Fagron	Urea	9K27-M01	1
Fagron	Urea	16D22-B02	3
Fagron	Urea	16D22-B02-322532	4
Fagron	Urea	16D22-B02-3221532	1
Fagron	Urea	2016142709	4
Fagron	Urea	18-B03-303581	1
Fagron	Urea	16D22-B02-322426	1
Fagron	Urea	16D22-B02-322526	4
Fagron	Uridine	11G07-N04	1
Fagron	Uridine	22041301	1
Fagron	Uridine	13D29-N06	3
Fagron	Uridine	16101301	1
Fagron	Uridine	13D08-N26	1
Fagron	Uridine	14F02-B09-305467	7
Fagron	Uridine	14F02-B09	1
Fagron	Uridine	721110090	4
Fagron	Uridine	12091402	1
Fagron	Uridine	14051405	1
Fagron	Uridine	13K18-B31	4
Fagron	Valine	121259	1
Fagron	Valine	13B13-N05	1
Fagron	Vanillin	14G02-B99-296988	1
Fagron	Vanillin	8L05-N07	1
Fagron	Vanillin	11G28-N09	1
Fagron	Vanillin	314Q2548	1
Fagron	Vanillin	232050-0004	1
Fagron	Vanillin	13J01-N07	6
Fagron	Wax, bleached	1701947	1
Fagron	Wax, bleached	15292022	7
Fagron	Wax, bleached	16A28-B10-322948	1
Fagron	Wax, bleached	2155	1
Fagron	Wax, yellow	13L02-B01-291452	1
Fagron	Wax, yellow	13B04-N04	2
Fagron	Wax, yellow	12D23-N05	2
Fagron	Wax, yellow	13L02-B01-294114	1
Fagron	Wax, yellow	13L02-B01-294113	2
Fagron	Wax, yellow	13102-B01	1
Fagron	Wild pansy, cut		1
Fagron	Wool alcohols	16235702	1
Fagron	Xylitol	13H28N10	1
Fagron	Xylitol	12J25-N06	2
Fagron	Xylitol	16c14-B02-320887	1
Fagron	Xylitol	14G04-B04-303524	2
Fagron	Xylitol	13/23-N06	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fagron	Xylitol	13H28-N10	2
Fagron	Xylometazoline hydroc...	1318-B03-295479	1
Fagron	Xylometazoline hydroc...	10416XMHC	1
Fagron	Xylometazoline hydroc...	12G31-N01	1
Fagron	Xylometazoline hydroc...	13I09-N02	1
Fagron	Xylometazoline hydroc...	13E06-N11	1
Fagron	Xylose	10121401	1
Fagron	Xylose	28111405	1
Fagron	Xylose	13133415	1
Fagron	Xylose	80414PJK01	1
Fagron	Xylose	14E-B03_300658	1
Fagron	Yarrow herb, cut	3613M-0165	1
Fagron	Yin Chen Hao	G034H1030421	2
Fagron	Zinc gluconate	14J30-B06-300520	1
Fagron	Captopril	14B26-B01-310127	6
Fagron	Mannitol	17C28-B05-33361	1
falcento	D-galactose	322619PF5	1
Falcento	Chromium (H1012)	16000520	5
Falcento	D-(+)-mannose	R-4638	1
Falcento	D-(+)-mannose	R-4639	1
Falcento	D-(+)-mannose	R-4641	1
Falcento	D-(+)-mannose	R4640	1
Falcento	D-galactose	2223PF4	3
Falcento	D-galactose	82619PF64743/16	1
Falcento	D-galactose	10440003	1
Falcento	D-galactose	12223PF4	2
Falcento	D-galactose	233619PF5	1
Falcento	D-galactose	5223PF4	3
Falcento	D-galactose	222PF4	1
Falcento	D-galactose	2014	1
Falcento	D-galactose	92619PF64744/16	1
Falcento	D-galactose	1223PF4	3
Falcento	D-galactose	13223PF4	1
Falcento	D-galactose	82619PF7	1
Falcento	D-galactose	122PF4	1
Falcento	D-galactose	72619PF7	1
Falcento	D-ribose	JXCZSJ14-D-010	1
Falcento	D-ribose	201209702	1
Falcento	D-ribose	201206712	2
Falcento	Xylose	2013051801	1
Fargo	Prednisolone acetate	12025-N05	1
fargon	Atropine sulfate	13K01-N01	1
fargon	Estriol	12a26-n02	1
fargon	Magnesium sulfate hep...	12E22-N02	4
fargon	Sodium citrate	100772-114-91	1
fargon	Triamcinolone acetoneide	13C13-N20	15
Fargon	Arginine	121429	1
Fargon	Ascorbyl palmitate	19D25-B04-361254	1
Fargon	Betamethasone dipropi...	16C04-B033-24513	1
Fargon	Betamethasone valerate	13k29-b32-289364	1
Fargon	Betamethasone valerate	7890	1
Fargon	Bifonazole	14k10-b01-300541	1
Fargon	Bismuth gallath, basic	13/18-B13-294108	1
Fargon	Caffeine citrate	12J05-N06	4
Fargon	Ciclosporin	13D26-No7	1
Fargon	Citric acid	12B01-NN01	1
Fargon	Clindamycin hydrochlo...	11L20-No5	1
Fargon	Fructose	15I21-B01-312701	1
Fargon	Hard fat in pastilles	41319264	1
Fargon	L-leucine	8160410903	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fargon	L-threonine	2114720407	1
Fargon	Metronidazole	14c14-b02-291238	1
Fargon	Metronidazole	13D08-N22	7
Fargon	Mometasone furoate, m...	14h06-b07-298481	1
Fargon	Mometasone furoate, m...	13H13-N04	12
Fargon	Naphazoline hydrochlo...	13/03-B03-295105	2
Fargon	Prednisolone, microni...	FP01409P	1
Fargon	Prilocaine hydrochlor...	17H02-B05	1
Fargon	Propranolol hydrochlo...	M140808	1
Fargon	Propyphenazone	11H31-N07	1
Fargon	Saccharin sodium	12A13-N03	1
Fargon	Triamcinolone acetonide	Nu2526M	1
Fargon	Triclosan	14F23-B02-295948	6
Fargon	Triclosan	4394403	1
Fargon	Triclosan	12C21-N05	1
Fargon	Urea	14e07-b05-297165	1
Fargon / Gehe	Cinnamon, cut	15104201	1
Fargon GmbH & Co.KG	Lidocaine hydrochlori...	16C09-B03-322144	7
Fargon GmbH & Co.KG	Melatonin	16G06-B01-324743	2
Fargon Services NL	Tetracaine hydrochlor...	15K04-B03-322151	13
Fargon/17.07.2013/Now...	Triamcinolone acetoni...	12 07-N04	1
Fargon/AHD	Prednisolone, microni...	18A29-B02-348044	3
Fargon/Alliance HC	Betamethasone valerate	16B18-B02-322111	2
Fargon/Alliance HC	Sodium citrate	13/24-B13-287921	2
Fargon/Alliance HC	Triamcinolone acetonide	15B02-B03-303830	7
Fargon/AllianceH.	Betamethasone dipropi...	17H16B04-340537	1
Fargon/AllianceH.	Cholesterol	16G05-B03-328689	1
Fargon/Anzag	Triamcinolone acetonide	13/26-N11	10
Fargon/Apo.am Park	Naphazoline hydrochlo...	539580002	1
Fargon/Gehe	Bacitracin	C2280092	2
Fargon/Gehe	Isosorbide dinitrate ...	17C29-B01-339447	7
Fargon/Kehr	Alfatradiol	16G05-B01-332379	2
Fargon/Kehr	Budesonide, micronized	16J04-B01-326583	1
Fargon/Kehr	Fluocinolonacetamid, ...	14L-19-B04	1
Fargon/Kehr	Testosterone propionate	17B09-B04-335269	1
Fargon/Kehr	Trometamol	17A16-B04-337589	1
Fargon/Noweda	Fluocinolonacetamid, ...	13L13-B35	1
Fargon/Phoenix	Caffeine citrate	15A05-B03-302746	12
Fargon/Phönix	Bacitracin	15C12-B06-312975	1
Fargon/Sana	Benzocaine	11K08-290049	1
fargon/Sanacorp	Sodium sulfate decahy...	12H31-N04	1
Fargon/Sanacorp	Citric acid	30.06.2016	1
Fargon/Sanacorp	Eosin Y	13/11-No3	1
fargon	Vanillin	13/01-N07	1
Fargron	Benzocaine	13L06-B02-294112	5
Fargron	Diclofenac sodium	14G17-B08-297032	7
Fargron	Tretinoin	1402007-01	2
Fargron/AHD	Gentamicin sulfate	12G04-N16	1
Fargron/Noweda	Sweet fennel	2870A150810	1
Farmalabor	Acetylcysteine	P1701913-003	1
Farmalabor	Folic acid hydrate	P1701610	1
Farmalabor	L-glutathione	P1700811	1
Farmalabor	Melatonin	R1804925	1
Farmalabor	Melatonin	P1800033-000	1
Farmalabor	Pyridoxine hydrochlor...	P1701531-000	1
Farmalabor / Farmalabor	Calcium ascorbate	P1104931-000A	1
Farmalabor / Farmalabor	Melatonin	P1201556-000	1
Farmalabor / Farmalabor	Salicylic acid	P1404195-002	1
Faron	Aluminium chloride he...	13I13-B03-287298	1
Faron	Tetracaine base	15D01-B03-306258	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fgaron	Clobetasol propionate	15H26-B03	1
Fgron	Caffeine citrate	2013Q-01282	1
Fiebig	Betamethasone, micron...	12J01-120301-m1	1
Fiebig	Chloramphenicol	18017810	3
Fiebig	Diltiazem hydrochloride	15230514	7
Fiebig	Diltiazem hydrochloride	17043227	1
Fiebig	Diphenylcyclopropenone	17B13-B02-334769	1
Fiebig	Estriol	17133711	2
Fiebig	Heparin sodium	16190116C0	1
Fiebig	Hydrochlorothiazide	115231003	1
Fiebig	Hydrocortisone acetate	1605410	2
Fiebig	Linseed	16011405	3
Fiebig	Sorbic acid	168009	1
Fiebig 05.11.12 ek:9...	Bilberry fruit, dried	2910A120410	1
Fiebig 16.01.13 ek: 4...	Wheat starch	11350911	3
Fiebig 16.11.12 ek:2,35	Walnut bark, cut	711B111020	5
Fiebig 27.06.13 EK3,34	Sodium benzoate	1304024-01	5
Fiebig 27.06.13 EK47,81	Betamethasone valerate	B9/131	1
Fiebig ek: 36.06EUR 31...	Minoxidil	111202403	1
Fiebig ek:12.27 01.08...	Diltiazem hydrochloride	11/01-N08	2
Fiebig ek:9.84EUR 24.09...	Peppermint leaf, cut	2181A110512	5
Fiebig/ Dextropur	Glucose monohydrate	L7045	1
Fiebig/Calo	Tetracaine hydrochlor...	17256704	1
Finteler Apotheke	Clobetasol propionate	16K24-B04-333121	3
Finteler Apotheke	Macrogols 4000	17035607	5
Finteler Apotheke	Triamcinolone acetone...	17057910	2
Finteler Apotheke	Vanillin	17117201	3
Finze	Calcium gluconate	917857359	1
Finze	Centaury, cut	11097610	1
Finze	Cetyl alcohol	10321703	1
Finze	Diphenhydramine hydro...	9F05-N04	2
Finze	Elder flower, rubbed	9297415	1
Finze	Equisetum stem, cut	11076305	8
Finze	Estradiol benzoate	1003040-01	2
Finze	Ethacridine lactate m...	90601	1
Finze	Frangula bark, cut	249091	1
Finze	Knotgrass, cut	10237202	1
Finze	Linseed	3607E-02240	1
Finze	Lycopodium herb, cut	9378803	1
Finze	Male speedwell wort, ...	4661B080915	1
Finze	Matricaria flower	12361703	6
Finze	Matricaria flower	1290B101122	1
Finze	Menthol	9221710	1
Finze	Sage leaf, cut	2301A110204	7
Finze	Sodium hydrogen carbo...	Pr.Nr.4311Q-01892	1
Finze	Sorbic acid	9280207	2
Finze	St. John's wort, cut	9297303	1
Finze	Tetracaine hydrochlor...	2906009-02	1
Finze	Thiamine chloride hyd...	9216114	1
Finze	Urea	1264	1
Finze	Wax, yellow	10145002	1
Finze	Wild thyme, cut	11031204	1
fischar	Sodium hydrogen carbo...	6522101	6
Fischar	Citric acid	6493071	2
Fischar	Citric acid	6498011	2
Fischar	Citric acid	649061	2
Fischar	Tartaric acid	704061	2
Fischar/Gehe	Estriol	14H29-B02-298072	3
Fischar/Gehe	Prednisolone, microni...	14A30-B02-294968	1
Fischar/Phönix	Tartaric acid	704021	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Fischer/Gehe	Carbomers 35000 / 50000	1410001-02	2
Fischer/Noweda	Citric acid	6496091	2
Fischer/VDL	Sodium hydrogen carbo...	6521111	5
Frago/Gehe	Estradiol valerate	15C10-B0508561040	1
Frago/Gehe	Isosorbide dinitrate ...	17C29-B01-335070	10
Frago/Noweda	Diphenylcyclopropenone	16F02-B10-324263	2
Frago/Noweda	Estradiol valerate	16F09-B01-329621	1
Frago/Noweda	Estradiol valerate	16F09-B01	2
Frago/Noweda	Lidocaine hydrochlori...	6108-05-0	1
Frago/Noweda	Metronidazole	16I19-B04-329163	1
fragon	Beclometasone dipropi...	700721-0004	1
fragon	Fructose	FR1413072	2
Fragon	Betamethasone dipropi...	18B06-B05-352731	1
Fragon	Estradiol valerate	15C10-B05-5314840	1
Fragon	Lactose monohydrate	10706040	1
Fragon	Miconazole nitrate	18C23-B02-348857	1
Fragon	Ofloxacin	15H19-B01-343412	1
Fragon	Prednisolone, microni...	18C19-B07-351225	2
Fragon	Progesterone, microni...	13A02-N02	1
fragon/AHD	Propranolol hydrochlo...	13/13-N04	7
Fragon/Phoenix	Citric acid	14G04-B09-293976	1
Fragon/Phoenix	Hypromellose 4000	18D04-b03-349970	1
Fragon/Phoenix	Neomycin sulfate	15J09-B04-322943	1
Fragon/Phönix	Citric acid	18F25-B02-351463	1
FRAGON/Poenix	Hypromellose 4000	16B16-B07-18276	1
Fragon/Sanacorp	Bismuth gallath, basic	16K17-B04-335593	1
Fragon/Sanacorp	Bismuth gallath, basic	17/26-B01-340552	1
Fragon/Sanacorp	Bismuth gallath, basic	18A-29-B06-346303	1
Fragon/Sanacorp	Hydrocortisone	17E10-B04-335056	2
Fragon/Sanacorp	Mometasone furoate, m...	17D19-B02-339277	7
Fragon/Sanacorp	Mometasone furoate, m...	17D19-B02-342642	1
Fragon/Sanacorp	Tartaric acid	13G03-N14	5
FragonFiebig	Mannitol	18A30-B07-348872	2
fragon	Emulsifying cetostear...	16L13-b01-326511	1
Fragon	Prednicarbate, micron...	NV0485M	1
Fragon	Prednisolone acetate	K07A20131101M	1
Fragon/noweda	Clotrimazole	CLT/1601/0005B	1
Fragon/noweda	Clotrimazole	17A19-B04	1
Fragon/Noweda	Citric acid	15G16-B01-330449	8
Fragon/Noweda	Mannitol	17A12-B04-330729	1
Fresenius	Benzoic acid	80715102	1
Fresenius	Benzydamine hydrochlo...	30062015099	1
Fresenius	Octenidine dihydrochl...	60715101	1
Fresenius	Potassium sorbate	220214221	1
Fresenius	Wheat starch	13305103	1
Frey*Lau GmbH	Menthol	1305010-01	11
Frey+Lau	Menthol	926636	1
Frey+Lau GmbH	Menthol	1509045-02	10
Frey+Lau GmbH/Phönix	Thymol	1312022-01	3
Fritz Schiele	Isoleucine (L-) (H1045)	1805027701	2
Fujian	Gentamicin sulfate	FG1706157	1
Fujian Fukang Pharmac...	Gentamicin sulfate	12G04-N16	1
Fujian Fukang Pharmac...	Gentamicin sulfate	14E08-B03-294848	2
Gahlke	Aniseed, pounded	21401	1
Gahlke	Bitter fennel	20362	1
Gahlke	Frankincense	21949	1
Gahlke	Parsley root, cut	20872	1
Gahlke	Sorrel herb, cut	21479	1
Gahlke	Sweet woodruff herb, ...	22100	1
galke	Aniseed	19654	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
galke	Bitter fennel	19938	3
galke	Melissa leaf, cut	20867	1
Galke	Acorus root, peeled	18944	1
Galke	Agrimony, cut	2002	1
Galke	Alchemilla, cut	21141	1
Galke	Alchemilla, cut	19835	1
Galke	Alchemilla, cut	22291	1
Galke	Aniseed	23464	1
Galke	Aniseed	21278	1
Galke	Aniseed	22274	1
Galke	Asparagus root, cut	20178	1
Galke	Bean shells, cut	21219	1
Galke	Bearberry leaf, cut	19936	3
Galke	Birch leaf, cut	20904	1
Galke	Birch leaf, cut	21226	1
Galke	Bitter fennel	22176	1
Galke	Blackcurrants	21348	1
Galke	Blue bottle flower wi...	23131	2
Galke	Boldo leaf, cut	19855	1
Galke	Burdock root, cut	19802	1
Galke	Camphor	23624	1
Galke	Caraway fruit	21257	1
Galke	Caraway fruit	27450	1
Galke	Caraway fruit	23627	1
Galke	Caraway fruit, pounded	18943	1
Galke	Centaury, cut	9721	1
Galke	Chelidonium majus her...	20937	3
Galke	Common broom herb, cut	12358803	4
Galke	Common wormwood herb,...	21637	1
Galke	Couch grass rhizome, ...	18781	1
Galke	Dandelion herb with r...	20531	6
Galke	Dandelion herb with r...	22093	2
Galke	Dandelion root, cut	126302	1
Galke	Dandelion root, cut	20658	1
Galke	Devil's claw root	20282	2
Galke	Devil's claw root	21018	2
Galke	Dill fruit	22422	4
Galke	Elder flower, rubbed	21934	1
Galke	Elecampane rhizome, cut	21569	2
Galke	Equisetum stem, cut	19663	5
Galke	Equisetum stem, cut	22061	1
Galke	Euphrasy herb, cut	20262	2
Galke	Gentian root, cut	19192	2
Galke	Geum urbanum herb, cut	20881	1
Galke	Ginger, peeled	21766	2
Galke	Ginger, peeled	20633	1
Galke	Goat's rue, cut	8378	2
Galke	Hawthorn leaf with fl...	21126	3
Galke	Hay flower	21926	1
Galke	Hop strobile	23363	1
Galke	Hop strobile	20928	3
Galke	Hop strobile	20928	2
Galke	Horse chestnut flower...	2614	1
Galke	Hyssop herb, cut	20939	1
Galke	Iceland moss, cut	12369	1
Galke	Ladys bedstraw herb,...	21942	2
Galke	Lavender flower	21416	2
Galke	Lemon peel, cut	22975	1
Galke	Liquorice root, peeled	23538	1
Galke	Lungwort herb, cut	19511	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Galke	Lycopodium herb, cut	16913	3
Galke	Mallow flower, cut	21119	1
Galke	Mallow leaf, cut	21445	1
Galke	Mallow leaf, cut	17133	1
Galke	Mallow leaf, cut	21454	1
Galke	Manna, powdered	19784	1
Galke	Manna, powdered	20669	1
Galke	Marian thistle fruit	22187	4
Galke	Marjoram herb, rubbed	23241	1
Galke	Marjoram herb, rubbed	22401	1
Galke	Masterwort root, cut	23074	1
Galke	Masterwort root, cut	18411	1
Galke	Maté leaf, cut	22365	1
Galke	Matricaria flower	18471	5
Galke	Matricaria flower	21652	1
Galke	Matricaria flower	18803	1
Galke	Meadowsweet herb, cut	20358	1
Galke	Mistletoe, cut	20328	1
Galke	Mistletoe, cut	22012	2
Galke	Nettle leaf, cut	21238	1
Galke	Nettle leaf, cut	20956	1
Galke	Nettle leaf, cut	23302	1
Galke	Nutmeg, powdered (spi. ...)	22739	1
Galke	Oak bark, cut	17824	2
Galke	Orris root	21303	1
Galke	Peppermint leaf, cut	20975	3
Galke	Peppermint leaf, cut	22218	1
Galke	Peppermint leaf, cut	22783	3
Galke	Primrose flower	22612	2
Galke	Psyllium seed	20968	2
Galke	Psyllium seed	22989	7
Galke	Raspberry leaf, cut	18060	2
Galke	Raspberry leaf, cut	23532	1
Galke	Raspberry leaf, cut	19499	2
Galke	Raspberry leaf, cut	21244	1
Galke	Red clover flower	23128	2
Galke	Red-rose petal	22641	2
Galke	Restharrow root, cut	12700	1
Galke	Ribwort plantain, cut	20259	1
Galke	Ribwort plantain, cut	21598	1
Galke	Rosemary leaf, cut	21390	1
Galke	Rosemary leaf, cut	19813	1
Galke	Rosemary leaf, cut	22777	1
Galke	Sage leaf, cut	21015	3
Galke	Sanicle herb, cut	23097	1
Galke	Senna leaf, cut	21778	1
Galke	Shepherd's purse, cut	21719	2
Galke	Shepherd's purse, cut	19483	1
Galke	Silverweed, cut	12891	2
Galke	St. John's wort, cut	21556	1
Galke	St. John's wort, cut	21034	1
Galke	St. John's wort, cut	21990	4
Galke	St. John's wort, cut	23054	1
Galke	Tall golden rod, cut	21079	2
Galke	Thyme, rubbed	1833	1
Galke	Thyme, rubbed	20719	2
Galke	Thyme, rubbed	28719	1
Galke	Thyme, rubbed	19386	2
Galke	Tormentil	21091	1
Galke	Tormentil	21744	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Galke	Turmeric root, cut	23123	1
Galke	Turmeric root, cut	20184	1
Galke	Turmeric, javanese, cut	15540	1
Galke	Turmeric, javanese, cut	15540	1
Galke	Valerian root, cut	21798	4
Galke	Valerian root, cut	22070	1
Galke	Walnut leaf, cut	21061	1
Galke	White deadnettle flower	22704	1
Galke	White deadnettle flower	21691	4
Galke	White deadnettle herb...	22253	2
Galke	White deadnettle herb...	19706	1
Galke	Willow bark, cut	20338	1
Galke	Willowherb, small flo...	22347	4
Galke	Willowherb, small flo...	47502	1
Galke	Wormwood herb, cut	21330	1
Galke	Yarrow herb, cut	21295	2
Galke	Yarrow herb, cut	23735	1
Galke	Yarrow herb, cut	23318	1
Galke	Yarrow herb, cut	22048	7
galke,	Calendula flower	12028	1
galke,	Lime flower, cut	2098	1
Gall Pharma	Prasterone	16/02-B02	1
Gall Pharma / Phoenix	Taurine	179.66	1
Gall-Pharma GmbH	Citrulline	n.v.	1
Gall-Pharma GmbH/Noweda	Citrulline	242.20	1
gasser	Angelica root	950049	1
Gatt Koller	Fructose	1544/05141914	1
Gatt-Koller	Alfatradiol	1803210	2
Gatt-Koller	Alfatradiol	180077661	1
Gatt-Koller	Ambroxol hydrochloride	180600901	1
Gatt-Koller	Ascorbic acid	1802406	2
Gatt-Koller	Ascorbic acid	18020406	1
Gatt-Koller	Betaine hydrochloride	2970/07102910	1
Gatt-Koller	Betamethasone valerate	1800117	1
Gatt-Koller	Biotin	1634211	1
Gatt-Koller	Boric acid	1803542	1
Gatt-Koller	Caffeine	1803747	1
Gatt-Koller	Calcium ascorbate	1900123	1
Gatt-Koller	Calcium citrate tetra...	1802139	3
Gatt-Koller	Camphor	1800202	3
Gatt-Koller	Carbomers 35000 / 50000	1800201	1
Gatt-Koller	Cetostearyl alcohol	1803965	1
Gatt-Koller	Cetyl alcohol	1690/06162416	1
Gatt-Koller	Clioquinol	1802295	2
Gatt-Koller	Clotrimazole	1801870	1
Gatt-Koller	Clotrimazole	201714s9	1
Gatt-Koller	Cocoa butter	1803838	1
Gatt-Koller	Dequalinium chloride	1800212	1
Gatt-Koller	Dexamethasone	1801087	1
Gatt-Koller	Diphenhydramine hydro...	606030315	1
Gatt-Koller	Erythromycin	1803029	2
Gatt-Koller	Estradiol hemihydrate	1801604	2
Gatt-Koller	Heparin sodium	1901027	1
Gatt-Koller	Hydrocortisone	1804169	1
Gatt-Koller	Hydrocortisone	20172187	1
Gatt-Koller	Hydrocortisone acetate	1802876	1
Gatt-Koller	Hydrocortisone acetate	1734/05154415	1
Gatt-Koller	Hydroquinone	1803302	4
Gatt-Koller	Hydroxyethylcellulose...	1802575	1
Gatt-Koller	Ketamine hydrochloride	1803056	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gatt-Koller	Lactose monohydrate	1702456	1
Gatt-Koller	Lidocaine	1900471	1
Gatt-Koller	Lidocaine hydrochlori...	1801084	1
Gatt-Koller	Macrogols 4000	1802391	1
Gatt-Koller	Macrogols 4000	2792/07153315	1
Gatt-Koller	Macrogols 4000	1420/04153415	1
Gatt-Koller	Macrogols 4000	1804126	1
Gatt-Koller	Macrogols 4000	4643/01170417	2
Gatt-Koller	Magnesium hydrogen ci...	1802721	1
Gatt-Koller	Magnesium hydrogen ci...	18.03714	1
Gatt-Koller	Magnesium hydrogen ci...	1800208	1
Gatt-Koller	Magnesium hydrogen ci...	1803714	2
Gatt-Koller	Mannitol	1803827	2
Gatt-Koller	Methyl-4-hydroxybenzo...	1803308	1
Gatt-Koller	Metronidazole	497/03151515	1
Gatt-Koller	Metronidazole	606293914	1
Gatt-Koller	Metronidazole	497/03151015	1
Gatt-Koller	Metronidazole	497/03154015	1
Gatt-Koller	Metronidazole	3728/11162017	1
Gatt-Koller	Minoxidil	1803103	2
Gatt-Koller	Naphazoline hydrochlo...	1803834	1
Gatt-Koller	Nicotinamide	4456/10104110	1
Gatt-Koller	Nicotinamide	221406112411	1
Gatt-Koller	Nystatin	1900206	2
Gatt-Koller	Phenylephrine hydroch...	330010164116	1
Gatt-Koller	Polymyxin B sulfate	17H23B01344840	1
Gatt-Koller	Potassium citrate	1800758	1
Gatt-Koller	Propyl-4-hydroxybenzo...	1800302	1
Gatt-Koller	Propyphenazone	1800011	1
Gatt-Koller	Quinine hydrochloride	3500/08103810	1
Gatt-Koller	Quinine sulfate dihyd...	1801042	1
Gatt-Koller	Quinine sulfate dihyd...	1803149	4
Gatt-Koller	Riboflavin	1803007	1
Gatt-Koller	Salicylic acid	1803840	2
Gatt-Koller	Sodium ascorbate	1802317	3
Gatt-Koller	Sodium citrate	1073/04121412	1
Gatt-Koller	Testosterone propionate	606213811	1
Gatt-Koller	Testosterone propionate	1803006	1
Gatt-Koller	Testosterone propionate	3478/09153915	1
Gatt-Koller	Tetracaine hydrochlor...	1803034	2
Gatt-Koller	Tetracaine hydrochlor...	KNC1825J18	1
Gatt-Koller	Tetracaine hydrochlor...	1900461	6
Gatt-Koller	Tetracycline hydrochl...	1700729	1
Gatt-Koller	Tretinoin	606210909	1
Gatt-Koller	Tretinoin	1702180	1
Gatt-Koller	Tretinoin	1804193	1
Gatt-Koller	Trometamol	1801776	1
Gatt-Koller	Trometamol	1804167	1
Gatt-Koller	Urea	1800697	1
Gatt-Koller	Urea	1803268	3
Gatt-Koller	Vanillin	201770800	1
Gatt-Koller / Phoenix	Benzoic acid	2992/07131314	2
Gatt-Koller / Phoenix	Caffeine	5898/12135213	1
Gatt-Koller / Phoenix	Caffeine	5898/12135213	1
Gatt-Koller / Phoenix	Calcium gluconate	550/03141214	1
Gatt-Koller / Phoenix	Chloramphenicol	5615/12115211	1
Gatt-Koller / Phoenix	Ephedrine hydrochloride	18502141114	1
Gatt-Koller / Phoenix	Fructose	3156/0813	1
Gatt-Koller / Phoenix	Hypromellose 4000	14C05-B03-291647	2
Gatt-Koller / Phoenix	Lactose monohydrate	5580/12120813	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gatt-Koller / Phoenix	Procaine hydrochloride	4428/10122714	2
Gatt-Koller / Phoenix	Propyphenazone	419/02123112	1
Gatt-Koller / Phoenix	Quinine sulfate dihyd...	4928/11132214	1
Gatt-Koller / Phoenix	Quinine sulfate dihyd...	2503/07143314	1
Gatt-Koller / Phoenix	Sodium tetraborate	6187/01142114	1
Gatt-Koller / Phoenix	Sorbic acid	2575/07142914	1
Gatt-Koller / Phoenix	Sorbitol	1718/0513	1
Gatt-Koller / Phoenix	Sorbitol	1718/05134013	1
Gatt-Koller / Phoenix	Tannin	333508134613	2
Gatt-Koller / Phoenix	Tannin	3335/08134613	1
Gatt-Koller / Phoenix	Tartaric acid	3334/08131514	1
Gatt-Koller / Phoenix	Tetracaine hydrochlor...	80/02140614	2
Gatt-Koller/Kwizda	Calcium glycerophosph...	4007/09123812	1
GATT-KOLLER/PHOENIX	Ephedrine hydrochloride	915/03141314	2
GATT-KOLLER/PHOENIX	Glucose monohydrate	243/02141714	1
GATT-KOLLER/PHOENIX	Glucose monohydrate	243/02142614	1
GATT-KOLLER/PHOENIX	Methyl-4-hydroxybenzo...	6112/01142814	1
GATT-KOLLER/PHOENIX	Papaverine hydrochlor...	2927/07131314	1
GATT-KOLLER/PHOENIX	Propyl-4-hydroxybenzo...	5080/11130214	1
GATT-KOLLER/PHOENIX	Sodium salicylate, po...	4853/11124512	1
gebepfarm	Sodium monohydrogenph...	20160909-2	1
gehe	Clobetasol propionate	APCLPM03616C	1
gehe	Hydrocortisone acetate	17219707	23
gehe	Pregnenolone	Apc07160801	1
gehe	Salicylic acid	16I15-B08-331442	1
gehe	Spironolactone	17E22-B11-335455	1
Gehe	Acetylcysteine	16139813	1
Gehe	Acetylsalicylic acid	10365302	1
Gehe	Acorus root, peeled	6701B130807	9
Gehe	Alchemilla, cut	1341102	1
Gehe	Alfatradiol	16087501	7
Gehe	Allantoin	14247603	1
Gehe	Aluminium chloride he...	130701102	1
Gehe	Aluminium chloride he...	12114902	3
Gehe	Amifampridine	13F07-N07	1
Gehe	Angelica root	285732	1
Gehe	Aniseed	10247203	1
Gehe	Aniseed, pounded	13149903	1
Gehe	Arnica flower	11259801	2
Gehe	Atropine sulfate	16326802	3
Gehe	Bacitracin	C2280065	3
Gehe	Benzocaine	111214	1
Gehe	Betamethasone valerate	14H28-B01-302863	5
Gehe	Betamethasone valerate	12003409	4
Gehe	Betamethasone valerate	1712E-01181	1
Gehe	Betamethasone valerate	14H25-B01-302863	1
Gehe	Betamethasone valerate	13091211	4
Gehe	Bifonazole	15152703	1
Gehe	Bitter fennel	2860A120412	24
Gehe	Blackberry leaf, cut	10136409	1
Gehe	Caffeine citrate	13B08-N02	3
Gehe	Calcium pantothenate	1004024-01	1
Gehe	Calendula flower, cut	260205	1
Gehe	Camphor	17134702	2
Gehe	Carbomers 35000 / 50000	13194505	3
Gehe	Cetostearyl alcohol	15272713	1
Gehe	Cetyl palmitate 15	8209646	1
Gehe	Chloramphenicol	1107002-02	2
Gehe	Chlortetracycline hyd...	16B10-B09	5
Gehe	Chlortetracycline hyd...	16B10-B09-322716	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Ciclopirox olamine	14E21-M04	1
Gehe	Citric acid	12043111	1
Gehe	Citric acid	11C24-N03	8
Gehe	Clindamycin hydrochlo...	13a22-n01	1
Gehe	Clindamycin hydrochlo...	700168-0002	1
Gehe	Clindamycin hydrochlo...	11K04-N12	2
Gehe	Clindamycin hydrochlo...	11E12-N02	1
Gehe	Clobetasol propionate	10D21-N04	1
Gehe	Clobetasol propionate	12342310	11
Gehe	Clobetasol propionate	12107606	1
Gehe	Clotrimazole	1108053-01	2
Gehe	Clotrimazole	14C04-B31-294952	2
Gehe	Coltsfoot leaf, cut	2021B121204	9
Gehe	Couch grass rhizome, ...	283061	1
Gehe	Couch grass rhizome, ...	12186704	1
Gehe	Dandelion herb with r...	13019803	1
Gehe	Dandelion herb with r...	6611B110104	1
Gehe	Deoxycholic acid	2012060253	1
Gehe	Devil's claw root	1114005	1
Gehe	Dexamethasone	1203032-02	4
Gehe	Dexamethasone	15054629	9
Gehe	Diclofenac sodium	15278502	3
Gehe	Diltiazem hydrochloride	17337026	9
Gehe	Diltiazem hydrochloride	12133009	9
Gehe	Dithranol	1311029-01	2
Gehe	Elder berries	12267807	1
Gehe	Elder flower, rubbed	1250213	1
Gehe	Elder flower, rubbed	1090413	1
Gehe	Elecampane rhizome, cut	11196906	1
Gehe	Elecampane rhizome, cut	12113302	2
Gehe	Emulsifying cetostear...	1701031-01	4
Gehe	Equisetum stem, cut	11076301	4
Gehe	Equisetum stem, cut	12288803	13
Gehe	Equisetum stem, cut	10375515	1
Gehe	Erythromycin	11D13-N09	3
Gehe	Erythromycin	16J24-B10-336735	7
Gehe	Erythromycin	18F05-B01-352501	1
Gehe	Estradiol benzoate	11L23-N01	2
Gehe	Estradiol hemihydrate	82404160	1
Gehe	Estradiol hemihydrate	15038016	1
Gehe	Fenugreek	11313204	1
Gehe	Fenugreek	13013811	2
Gehe	Frangula bark, cut	11254004	4
Gehe	Frangula bark, cut	11254003	2
Gehe	Frangula bark, cut	11254006	4
Gehe	Fumitory, cut	12103402	2
Gehe	Galeopsis herb, cut	3911A130320	6
Gehe	Gentamicin sulfate	15226517	2
Gehe	Gentamicin sulfate	FG1307194	1
Gehe	Glucose monohydrate	13157026	6
Gehe	Glycerol monostearate...	14240605	2
Gehe	Hard fat in pastilles	13243403	9
Gehe	Hawthorn leaf with fl...	1351B111207	2
Gehe	Horse chestnut leaf, ...	2551B131024	1
Gehe	Hydrochlorothiazide	12E04-N02	2
Gehe	Hydrocortisone	13H08-B01-285471	1
Gehe	Hydrocortisone	13069510	7
Gehe	Hydrocortisone	17E10-B04-335656	2
Gehe	Hydrocortisone	13H08-B01-285469	2
Gehe	Hydrocortisone acetate	14043905	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Hydroquinone	12156706	3
Gehe	Hydroxyethylcellulose...	12156908	5
Gehe	Hydroxyethylcellulose...	111002404	1
Gehe	Hydroxyethylcellulose...	13A18-N01	5
Gehe	Iceland moss, cut	11053915	1
Gehe	Iceland moss, cut	11101603	1
Gehe	Ivy leaf, cut	2061B080402	1
Gehe	Juniper, pounded	13020903	1
Gehe	Kaolin, heavy	11350107	3
Gehe	Ketoconazole	11070019	1
Gehe	Lidocaine	14E13-B05-295552	5
Gehe	Lidocaine hydrochlori...	12C05-N02	3
Gehe	Lidocaine hydrochlori...	1190313	1
Gehe	Lidocaine hydrochlori...	LDH/116023	2
Gehe	Liquorice root, peeled	270853	3
Gehe	Liquorice root, peeled	268835	1
Gehe	Liquorice root, peeled	1250213	1
Gehe	Liquorice root, natur...	277004	2
Gehe	Liquorice root, natur...	283316	1
Gehe	Liquorice root, natur...	284243	2
Gehe	Macrogols 1500	14072409	2
Gehe	Mallow flower, cut	1501B110207	1
Gehe	Mallow flower, cut	1501B120508	6
Gehe	Mallow leaf, cut	10290512	1
Gehe	Mallow leaf, cut	2131B110420	2
Gehe	Mannitol	17L04-B03-345230	2
Gehe	Mannitol	12E24-N08	2
Gehe	Mannitol	13E24-N01	1
Gehe	Mannitol	12J26-N15	2
Gehe	Marian thistle herb, ...	10365602	2
Gehe	Marshmallow root, pee...	6071A120306	2
Gehe	Matricaria flower	11214802	1
Gehe	Matricaria flower	11255204	1
Gehe	Matricaria flower	12337605	1
Gehe	Matricaria flower	12133009	1
Gehe	Melatonin	12G19-N03	4
Gehe	Melatonin	13E22-N12	2
Gehe	Melatonin	11/26-N07	1
Gehe	Melatonin	11L21-N04	1
Gehe	Melilot herb, cut	110324	1
Gehe	Melissa leaf, cut	12032314	4
Gehe	Methyl-4-hydroxybenzo...	11090020	1
Gehe	Metronidazole	10368318	2
Gehe	Metronidazole	16H15-B50-326734	2
Gehe	Minoxidil	15436406	5
Gehe	Minoxidil	AP006/13M	3
Gehe	Mometasone furoate, m...	16353201	7
Gehe	Motherwort herb, cut	12231502	3
Gehe	Neomycin sulfate	183360-0009	1
Gehe	Neomycin sulfate	194A-01926	1
Gehe	Neostigmine bromide	12172307	1
Gehe	Neostigmine bromide	12172301	2
Gehe	Nettle leaf, cut	2381A120619	13
Gehe	Nettle leaf, cut	10162701	1
Gehe	Nettle root, cut	6631B120927	5
Gehe	Nystatin	1207009-01	10
Gehe	Oxytetracycline hydro...	12g10-n04	1
Gehe	Parsley root, cut	6441A110210	3
Gehe	Peppermint leaf, cut	2181B130717A	1
Gehe	Peppermint leaf, cut	2181B130717B	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Pilocarpine hydrochlo...	15198610	2
Gehe	Polymyxin B sulfate	170101001	5
Gehe	Potassium hydrogen ta...	11345808	1
Gehe	Potassium sorbate	11261304	3
Gehe	Prednicarbate, micron...	12B14-N02	1
Gehe	Prednisolone acetate	13262609	3
Gehe	Prednisolone, microni...	12C22-N05	8
Gehe	Prednisolone, microni...	13A07-N14	8
Gehe	Prednisolone, microni...	12324216	12
Gehe	Prednisolone, microni...	13E15-No6	1
Gehe	Prednisolone, microni...	13343107	2
Gehe	Prednisolone, microni...	16C15-B01	2
Gehe	Prednisolone, microni...	16C-15-B02	1
Gehe	Prednisolone, microni...	16E09-B07	1
Gehe	Prednisolone, microni...	14A30-B02-294968	6
Gehe	Prednisone	12L14-M03	1
Gehe	Prednisone	12232714	2
Gehe	Primrose flower	11374401	1
Gehe	Propranolol hydrochlo...	16E09-B09	2
Gehe	Propranolol hydrochlo...	j051201.1	1
Gehe	Resocinol	11L13-N02	1
Gehe	Rhubarb, cut	13235208	2
Gehe	Ribwort plantain, cut	1080413	1
Gehe	Rue herb, cut	12148802	1
Gehe	Salicylic acid	1050313	1
Gehe	Salicylic acid	14L02-B02-303558	12
Gehe	Salicylic acid	11D15-N08	1
Gehe	Salicylic acid	165023-02	1
Gehe	Sildenafil citrate	13F24-N03	1
Gehe	Sildenafil citrate	16049805	3
Gehe	Silverweed, cut	11328101	1
Gehe	Sodium dodecyl sulfate	14127404	1
Gehe	Sodium hydrogen carbo...	12k09-N05	1
Gehe	Sodium hydrogen carbo...	11078315	1
Gehe	Sodium hydrogen carbo...	1203030-02	8
Gehe	Sodium monohydrogenph...	12109001	2
Gehe	Sodium salicylate, po...	13G31-B03-287961	1
Gehe	Sodium sulfate decahy...	12257902	6
Gehe	Spironolactone	12049105	3
Gehe	St. John's wort, cut	9297305	1
Gehe	St. John's wort, cut	11185602	1
Gehe	Sulfacetamide	1112032-01	7
Gehe	Sulfacetamide	11085904	1
Gehe	Sulfacetamide	1505025-01	1
Gehe	Tall golden rod, cut	11280311	2
Gehe	Tall golden rod, cut	11207011	1
Gehe	Tannin	12336016	9
Gehe	Tartaric acid	7284696	1
Gehe	Tartaric acid	1325-B01-285487	1
Gehe	Tartaric acid	703061	1
Gehe	Tartaric acid	16146806	3
Gehe	Tetracycline hydrochl...	1110029-01	3
Gehe	Thyme, rubbed	4622A110120	1
Gehe	Triamcinolone acetonide	15B02-B03-303827	3
Gehe	Triamcinolone acetonide	13I26-N11	1
Gehe	Triamcinolone acetonide	700124-0007	1
Gehe	Triamcinolone acetonide	12301510	6
Gehe	Triamcinolone acetonide	13124805	4
Gehe	Triclosan	12106-n01	3
Gehe	Triclosan	13391811	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe	Triclosan	13391804	12
Gehe	Trometamol	11160708	1
Gehe	Trometamol	11025103	1
Gehe	Tropicamide	15F18-B04	1
Gehe	Tropicamide	15F-18B-04	1
Gehe	Verbena herb, cut	4651B110310	2
Gehe	Walnut bark, cut	711B110128	3
Gehe	Walnut leaf, cut	2111B151028	2
Gehe	Walnut leaf, cut	10050412	1
Gehe	Wax, yellow	12079412	1
Gehe	Wild pansy, cut	11216103	3
Gehe	Wild thyme, cut	1090413	1
Gehe	Wormwood herb, cut	11285811	1
Gehe	Wormwood herb, cut	13381901	1
Gehe	Xylometazoline hydroc...	10118802	1
Gehe	Yarrow herb, cut	1128503	1
Gehe	Yarrow herb, cut	4181A111122	5
Gehe	Yarrow herb, cut	279023	3
GEHE	Arginine	17G10-B07-341856	2
GEHE	Betamethasone dipropi...	15412914	6
GEHE	Camphor	16235408	2
GEHE	Captopril	17210203	1
GEHE	Chloramphenicol	160200202	1
GEHE	Clotrimazole	177111613	1
GEHE	Estradiol valerate	16F09-B01-332757	2
GEHE	Glucose	18012902	11
GEHE	Hydrochlorothiazide	18017302	2
GEHE	Ibuprofen	15103011	1
GEHE	Iodoform	17086003	1
GEHE	Mannitol	17B08-B04-334312	1
GEHE	Neomycin sulfate	1502001-03	1
GEHE	Paraffin, hard	15295409	1
GEHE	Potassium sorbate	17170008	15
GEHE	Sildenafil citrate	17344003	3
GEHE	Sodium citrate	17A26-B04-332866	6
GEHE	Tetracaine base	16G29-B02-325569	4
GEHE	Triamcinolone acetone...	1710002-01	35
GEHE	Triamcinolone acetone...	16K07-B06-332103	1
Gehe / 08.03.13 Redwood	Blessed thistle, cut	205007	1
Gehe / 26.10.12 Fragon	Hypromellose 4000	9D30-N05	1
Gehe 01.09.2014	Sorbitol	13278217	1
GEHE 02.05.2013	Prilocaine hydrochlor...	12J26.N09	1
GEHE 02.05.2014	Tetracaine base	14A16-B04-290208	5
GEHE 02.12.2014	Lidocaine	13J30-B02-292220	2
GEHE 04.08.2014	Clioquinol	13114418	7
GEHE 06.05.2014	Macrogols 4000	13356608	3
Gehe 06.10.14	Carmellose sodium (40...	13299312	1
Gehe 07.12.18	Mometasone furoate, m...	1711027_02	1
gehe 07.12.2013	Sage leaf, cut	2301/7120403	1
GEHE 10.03.2014	Oak bark, cut	761A11816	1
gehe 10.12.13	Mannitol	1305043/01	1
gehe 10.12.13	Raspberry leaf, cut	229117121120	1
GEHE 21.10.2014	Deoxycholic acid	13F13-N06	1
GEHE 22.03.2014	Melatonin	13167302	1
Gehe 23.11.18	Metronidazole	1706022-02	2
GEHE 27.05.2014	Ketamine hydrochloride	13H19-N01	1
Gehe 27.11.14	Glucose	14142309	1
GEHE 28.11.2014	Sorbitol	13278245	1
Gehe am 03.02.2014	Betamethasone valerate	13091210	7
Gehe am 04.03.14	Ephedrine hydrochloride	13323303	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Gehe am 24.2.14	Gentamicin sulfate	1305047-01	8
Gehe am 27.03.2014	Juniper, pounded	2890A130513	1
Gehe Kassel	Calcium citrate tetra...	14133307	2
Gehe Kassel 26.06.2015	Propranolol hydrochlo...	15082401	6
Gehe/ 27.10.12 EuRho	Menthol	11219216	1
Gehe/Audar	Dithranol	APDT3752012D	2
Gehe/Euro	Sodium citrate	1412014-01	20
Gehe/Euro	Sodium citrate	1404033-01	9
Gehe/EuRo	Triamcinolone acetonide	1801020-01	27
Gehe/Fischer	Tartaric acid	706061	1
Gehe/Fragon	Prednisolone, microni...	16E09B07328150	1
Gehe/Fragon	Testosterone propionate	17B09B03332376	2
Gehe/Fragon	Triclosan	15C10B06312981	2
Gehe/Gehe	Ascorbic acid	L1505001	1
Gehe/Gehe	Hydroquinone	14280110	1
Gehe/Henry Lamotte	Tartaric acid	17/12	1
Gehe01/2018	Metronidazole	14C14-B02-291237	1
Gehe-Eurho 22.11.2017	Tretinoin	2171127Ma7	1
Gehe-Fragon	Prednisolone, microni...	14B18B05294966	1
Gehe-Fragon	Triclosan	15G07B03310036	2
gepepharm/Phönix	Hydroxyethylcellulose...	160603002	1
Gesundform	Ascorbic acid	L1110704	1
Gesundform	Mistletoe, cut	216482	1
ghfgh	Potassium citrate	3654645	1
Giusto Faravelli	D-(+)-mannose	161004/1710215001	1
gonmisol	Nicotinic acid	17-18/NCN(P)/B/008	1
Gutenberg-Apotheke	Agnus castus fruit	61765	1
Gutenberg-Apotheke	Aniseed	1333902	1
Gutenberg-Apotheke	Carmellose sodium (40...	4811E-03302	1
Gutenberg-Apotheke	Chloramphenicol	13A31-B02-286710	1
Gutenberg-Apotheke	Ground ivy herb, cut	11374610	3
Gutenberg-Apotheke	Hydroxyethylcellulose...	11048105	1
Gutenberg-Apotheke	Lidocaine	13J30-B02	3
Gutenberg-Apotheke	Macrogols 4000	13356609	2
Gutenberg-Apotheke	Mannitol	13189221	1
Gutenberg-Apotheke	Nettle leaf, cut	18897	1
Gutenberg-Apotheke	Oak bark, cut	11351206	13
Gutenberg-Apotheke	Raspberry leaf, cut	28911	1
Gutenberg-Apotheke	Rice starch	12360216	1
Gutenberg-Apotheke	Salicylic acid	12374208	3
Gutenberg-Apotheke	Urea		5
Haihang	Prasterone	20130507	1
Haihang Industry	Estriol	20130824	1
Haihang Industry	Prasterone	20130816	4
Haihang Industry	Pregnenolone	20140329	1
Haihang Industry	Testosterone	20140114-1	1
Haihang Industry	Testosterone	20140114-2	1
hannes	Taurine	1202125	1
Hanns G.Werner GmbH	Sucrose	8731158	1
Haushalt	Sucrose	123456	1
Haustee	Mistletoe, cut	57010507	1
HBK Gesundheit	Maltodextrin	1143	3
HBK Gesundheit	Maltodextrin	2343	1
HCK	Aminomix NAC (H1002)	16001453	1
HCK	Arginine (L-) (H1036)	14000564	1
HCK	Bioflavonoid complex ...	17000256	3
HCK	Calcium (H1008)	15001672	3
HCK	Carnosine (L-) (H1014)	1813979202	1
HCK	Clobetasol propionate	11/05-N06	1
HCK	Coenzyme Q10 (H1013)	16001598/1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
HCK	Creatine (H1031)	17001177	2
HCK	Ginseng (H1018)	14002047	1
HCK	Glycine (H1021)	903002.1803	1
HCK	Inositol hexanicotina...	19000750/4	1
HCK	Iron (H1015)	1259	1
HCK	Isoleucine (L-) (H1045)	15001641	3
HCK	L-5-HTP (H1035)	16000441	1
HCK	Leucine (L-) (H1046)	965101	1
HCK	Lycopene (H1058)	17000619/0	2
HCK	Methylsulfonylmethane...	1400921/1	1
HCK	Paraffin, hard	16277705	3
HCK	Silymarin (H1071)	15001046/0	1
HCK	Threonine (L-) (H1052)	14001718/9	1
HCK	Tryptophan (L-) (H1053)	15000441	1
HCK	Tryptophan (L-) (H1053)	16001404	1
HCK	Tyrosine (L-) (H1054)	1380	1
HCK	Vitamin B12 (H1102)	15000993/0	2
HCK	Vitamin C (H1104)	16000203/0	3
HCK	Vitamin E NAT (H1107)	14001822/0	1
hedinger	Hydroxyethylcellulose...	1606030	1
Hedinger	Bismuth gallath, basic	16038310	5
Hedinger	Captopril	14B26-B40	3
Hedinger	Dandelion herb with r...	10121614	1
Hedinger	Erythromycin	15245002	19
Hedinger	Linseed	16011402	1
Hedinger	Macrogols 4000	13G19-N08	1
Hedinger	Magnesium sulfate hep...	204036	1
Hedinger	Magnesium sulfate hep...	204037	2
Hedinger	Miconazole nitrate	16144720	7
Hedinger	Minoxidil	14B04-B07	3
Hedinger	Rhubarb, powdered	13193110	1
Hedinger	Shu Di (Huang)	g210h1601221	4
Hedinger	Sodium benzoate	13I27-B03	1
Hedinger	Threonine (L-) (H1052)	25708	6
Hedinger / Phönix	Citric acid monohydrate	202029	4
Hedinger/ Sanacorp	Clobetasol propionate	APAPL/035/E-12C	3
Hedinger/Alliance H.	Clotrimazole	15395617	9
Hedinger/Ebert	Vanillin	13A30-N10	1
Hedinger/Phönix	Sodium sulfate decahy...	205025	1
Heidelberger Chlorella	Magnesium hydrogen ci...	90196	1
Heidelberger Chlorell...	Magnesium hydrogen ci...	10032	1
Heinrich Kle/Sanacor	Linseed	7307A170111	2
Henry Lamotte	Macrogols 4000	1213test	1
Henry Lamotte Oils GmbH	Lidocaine	12 21-N06	1
Henry Lamotte Oils GMBH	Macrogols 4000	16267308/170127M	1
Hepart	Aminomix NAC (H1002)	18001998/0	1
Hepart	Aminomix NAC (H1002)	16001453/0	2
Hepart	Aminomix NAC (H1002)	17000766/0	2
Hepart	Antioxidants NAT comp...	18001182	1
Hepart	Antioxidants NAT comp...	17001754	1
Hepart	Antioxidants NAT comp...	17000265	1
Hepart	Antioxidants NAT comp...	16000968	1
Hepart	Arginine (L-) (H1036)	17000878/1	2
Hepart	Arginine (L-) (H1036)	17000878	4
Hepart	Arginine (L-) (H1036)	17000877/0	2
Hepart	Bioflavonoid complex ...	14001675/5	1
Hepart	Biotin (H1007)	14000862/4	2
Hepart	Calcium (H1008)	15001667/0	1
Hepart	Calcium (H1008)	17000204	1
Hepart	Calcium (H1008)	150016740	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart	Calcium (H1008)	170002040	1
Hepart	Calcium (H1008)	17000205	1
Hepart	Calcium (H1008)	17000206	1
Hepart	Calcium (H1008)	17000251/0	1
Hepart	Carnitine (L-) (H1038)	16000453	1
Hepart	Carnitine (L-) (H1038)	16000454/0	7
Hepart	Carnitine (L-) (H1038)	160004531	3
Hepart	Carnitine (L-) (H1038)	16000453/0	1
Hepart	Carnitine (L-) (H1038)	16000451/0	1
Hepart	Carnitine (L-) (H1038)	16000454	1
Hepart	Carnitine (L-) (H1038)	15000054/1	2
Hepart	Carnitine (L-) (H1038)	160000453	1
Hepart	Carnosine (L-) (H1014)	18001109/0	1
Hepart	Catechin extract (H10...	17000409/0	1
Hepart	Catechin extract (H10...	115000206/5	1
Hepart	Chondroitin sulfate ...	17000211/0	1
Hepart	Chondroitin sulfate ...	14001985/5	1
Hepart	Chromium (H1012)	16000520/3	3
Hepart	Chromium (H1012)	160005203	1
Hepart	Chromium (H1012)	16000520/8	5
Hepart	Coenzyme Q10 (H1013)	16000402/1	1
Hepart	Coenzyme Q10 (H1013)	16001599	2
Hepart	Copper (H1032)	150002117	1
Hepart	Copper (H1032)	15000211/5	1
Hepart	Curcumin and piperine...	17000258	2
Hepart	Curcumin and piperine...	17000259/0	1
Hepart	Curcumin and piperine...	1017	1
Hepart	Curcumin and piperine...	180016870	1
Hepart	Garcinia cambogia (H1...	16001252/1	1
Hepart	Ginkgo (H1024)	17001832	1
Hepart	Ginkgo (H1024)	18001110	1
Hepart	Ginkgo (H1024)	150010391	1
Hepart	Ginseng (H1018)	18000794/0	1
Hepart	Glucosamine sulfate (...)	170000073	1
Hepart	Glucosamine sulfate (...)	34001	1
Hepart	Glucosamine sulfate (...)	17001245	2
Hepart	Glutamine (L-) (H1043)	17000901	1
Hepart	Glutamine (L-) (H1043)	17000901/1	1
Hepart	Glutathione (H1020)	16000938/0	2
Hepart	Glutathione (H1020)	15000478	2
Hepart	Glycine (H1021)	17001660/0	1
Hepart	Iron (H1015)	15001258	2
Hepart	Iron (H1015)	15001259/2	2
Hepart	Isoflavones (H1025)	16001644	1
Hepart	Isoleucine (L-) (H1045)	17000900	1
Hepart	L-5-HTP (H1035)	18000913	1
Hepart	Leucine (L-) (H1046)	18001305	1
Hepart	Lipoic acid (α-) (H10...	18000250	1
Hepart	Lycopene (H1058)	17000619	2
Hepart	Lysine (L-) (H1047)	16001395	3
Hepart	Lysine (L-) (H1047)	16001503	1
Hepart	Magnesium (H1060)	16000269/0	1
Hepart	Magnesium (H1060)	16001704/0	3
Hepart	Magnesium (H1060)	16001704	1
Hepart	Methylsulfonylmethane...	17001247	1
Hepart	OPC Grape seed (H1066)	16001439	3
Hepart	OPC Grape seed (H1066)	17001805	2
Hepart	OPC Grape seed (H1066)	17001805/1	1
Hepart	Phenylalanine (L-) (H...	17000876	1
Hepart	Phenylalanine (L-) (H...	170008760	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart	Phenylalanine (L-) (H...	17000876/1	3
Hepart	Potassium (H1030)	16001645	1
Hepart	Potassium (H1030)	15000599/1	1
Hepart	Potassium (H1030)	16001646	1
Hepart	Proline (-L) NF (H1050)	17000952	2
Hepart	Selenium (H1067)	16001350	1
Hepart	Selenium (H1067)	160013497/0	1
Hepart	Selenium (H1067)	15000158/1	1
Hepart	Selenium yeast (H1068)	71901	2
Hepart	Taurine (H1084)	160009470	1
Hepart	Taurine (H1084)	16000947	4
Hepart	Taurine (H1084)	16000946	2
Hepart	Tetracaine base	15G30-B04-310590	1
Hepart	Threonine (L-) (H1052)	17000875	1
Hepart	Trace elements JK (H1...	15000529/0	1
Hepart	Trace elements JK (H1...	17000869	2
Hepart	Trace elements JK (H1...	17000869/0	2
Hepart	Tryptophan (L-) (H1053)	16001394/0	1
Hepart	Tyrosine (L-) (H1054)	18000412/0	2
Hepart	Vitamin B1 (H1101)	17000914	1
Hepart	Vitamin B6 (H1103)	15000879	1
Hepart	Vitamin B6 (H1103)	554938001	1
Hepart	Vitamin C (H1104)	170009150	1
Hepart	Vitamin C (H1104)	16000231/0	2
Hepart	Vitamin complex (H1100)	17001751/0	2
Hepart	Vitamin complex (H1100)	17000880	1
Hepart	Vitamin complex (H1100)	16000958	1
Hepart	Vitamin complex (H1100)	17001654	1
Hepart	Vitamin complex (H1100)	17001750/0	2
Hepart	Vitamin D3 (H1105)	17000874	2
Hepart	Vitamin D3 (H1105)	16001212/0	1
Hepart	Vitamin E NAT (H1107)	17000892/0	1
Hepart	Vitamin E NAT (H1107)	16000591/0	2
Hepart	Vitamin E NAT (H1107)	160005910	1
Hepart	Vitamin E NAT (H1107)	17000890	1
Hepart	Zinc (H1110)	16000519	1
Hepart	Zinc (H1110)	16000518	3
Hepart / Unisan	Calcium (H1008)	17000250	1
Hepart / Unisan	Catechin extract (H10...	170004092	1
Hepart / Unisan	Coenzyme Q10 (H1013)	160015991	1
Hepart / Unisan	Lipoic acid (α-) (H10...	150015262	1
Hepart / Unisan	Ornithine (L-) (H1048)	15000004	2
Hepart / Unisan	Potassium (H1030)	160016460	1
Hepart / Unisan	Vitamin C (H1104)	170009170	1
hepart ag	Silymarin (H1071)	1244950401	1
hepart ag	Vitamin complex (H1100)	1015954801	1
hepart ag	Vitamin complex (H1100)	10159548011	1
Hepart AG	Acacia gum (H1001)	728201	4
Hepart AG	Acacia gum (H1001)	1211728201	3
Hepart AG	Alchemilla, cut	21963	3
Hepart AG	Alchemilla, cut	B-698kbA	1
Hepart AG	Aminomix NAC (H1002)	1054019101	3
Hepart AG	Aminomix NAC (H1002)	15000320	2
Hepart AG	Aminomix NAC (H1002)	16001452	1
Hepart AG	Aminomix NAC (H1002)	16001552/0	1
Hepart AG	Aniseed	B216346,kbA	1
Hepart AG	Antioxidants NAT comp...	1034030902	1
Hepart AG	Antioxidants NAT comp...	16000376	3
Hepart AG	Antioxidants NAT comp...	1037030904	1
Hepart AG	Antioxidants NAT comp...	17001754/0	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart AG	Antioxidants NAT comp...	1034019901	1
Hepart AG	Antioxidants NAT comp...	15000878/0	2
Hepart AG	Antioxidants NAT comp...	52201	3
Hepart AG	Antioxidants NAT comp...	17000266/0	1
Hepart AG	Antioxidants NAT comp...	17000264/0	1
Hepart AG	Antioxidants NAT comp...	18001182/1	1
Hepart AG	Arginine (L-) (H1036)	14000563/2	1
Hepart AG	Arginine (L-) (H1036)	17000878/0	9
Hepart AG	Bioflavonoid complex ...	1284042501	2
Hepart AG	Bioflavonoid complex ...	42501	1
Hepart AG	Bioflavonoid complex ...	1284019001	2
Hepart AG	Bioflavonoid complex ...	68101	2
Hepart AG	Biotin (H1007)	29601	4
Hepart AG	Bitter fennel	B492, kbA	1
Hepart AG	Bitter orange flower ...	B120822	1
Hepart AG	Black cumin seed	20720, kbA	1
Hepart AG	Buckwheat herb, cut	B120509	2
Hepart AG	Calcium (H1008)	17000203/0	2
Hepart AG	Calcium (H1008)	17000205/0	2
Hepart AG	Calcium (H1008)	15001668/0	2
Hepart AG	Calcium (H1008)	14001933	1
Hepart AG	Calcium (H1008)	55501	8
Hepart AG	Calcium (H1008)	55401	6
Hepart AG	Calcium (H1008)	55201	1
Hepart AG	Calcium (H1008)	55001	3
Hepart AG	Caraway fruit	B-224, kbA	1
Hepart AG	Carnitine (L-) (H1038)	480g	1
Hepart AG	Carnitine (L-) (H1038)	16000453/1	2
Hepart AG	Carnitine (L-) (H1038)	26010	6
Hepart AG	Carnitine (L-) (H1038)	130605301	1
Hepart AG	Carnitine (L-) (H1038)	26001	1
Hepart AG	Carnitine (L-) (H1038)	1184026001	1
Hepart AG	Carnosine (L-) (H1014)	1812045901	2
Hepart AG	Carnosine (L-) (H1014)	1812979202	1
Hepart AG	Carnosine (L-) (H1014)	45901	2
Hepart AG	Catechin extract (H10...	15000206/2	2
Hepart AG	Catechin extract (H10...	1574018601	1
Hepart AG	Catechin extract (H10...	17000409/2	1
Hepart AG	Chondroitin sulfate ...	31801/1	1
Hepart AG	Chromium (H1012)	16000520/7	4
Hepart AG	Chromium (H1012)	14000429/4	1
Hepart AG	Chromium (H1012)	H101221	1
Hepart AG	Coenzyme Q10 (H1013)	130521-300	11
Hepart AG	Coenzyme Q10 (H1013)	1194015601	5
Hepart AG	Coenzyme Q10 (H1013)	16001599/1	2
Hepart AG	Coenzyme Q10 (H1013)	16001599/2	12
Hepart AG	Copper (H1032)	1544046401	1
Hepart AG	Copper (H1032)	13041501	7
Hepart AG	Copper (H1032)	1544007901	1
Hepart AG	Copper (H1032)	130415-K1	1
Hepart AG	Creatine (H1031)	15000749	1
Hepart AG	Curcumin and piperine...	13901	5
Hepart AG	Curcumin and piperine...	1534002201	1
Hepart AG	Curcumin and piperine...	15001589	1
Hepart AG	Daisy flower	A130307	1
Hepart AG	Dill fruit	L30057kbA	1
Hepart AG	Equisetum stem, cut	B-695, kbA	1
Hepart AG	Folic acid (H1016)	1613030201	6
Hepart AG	Folic acid (H1016)	15000125	1
Hepart AG	Folic acid (H1016)	52301	10

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart AG	Garcinia cambogia (H1...	7401	2
Hepart AG	Ginkgo (H1024)	15001038	1
Hepart AG	Ginkgo (H1024)	15001039/1	2
Hepart AG	Ginseng (H1018)	15001047	1
Hepart AG	Ginseng (H1018)	17000285/0	1
Hepart AG	Glucosamine sulfate (...)	1822038001	1
Hepart AG	Glucosamine sulfate (...)	38001	6
Hepart AG	Glutamine (L-) (H1043)	15000017/1	1
Hepart AG	Glutamine (L-) (H1043)	160009420	1
Hepart AG	Glutamine (L-) (H1043)	17000901/0	1
Hepart AG	Glutathione (H1020)	1254040201	1
Hepart AG	Glutathione (H1020)	16000938/2	1
Hepart AG	Glutathione (H1020)	1254025703	2
Hepart AG	Glycine (H1021)	17001660/1	1
Hepart AG	Glycine (H1021)	75501	1
Hepart AG	Glycine (H1021)	1803028101	1
Hepart AG	Glycine (H1021)	18000449/0	1
Hepart AG	Glycine (H1021)	14001730/2	1
Hepart AG	Hawthorn flower	13010401	1
Hepart AG	Inositol hexanicotina...	16000750/4	1
Hepart AG	Iron (H1015)	1114040901	4
Hepart AG	Isoflavones (H1025)	1314038801	1
Hepart AG	Isoflavones (H1025)	56701	3
Hepart AG	Isoflavones (H1025)	1314025704	1
Hepart AG	Isoleucine (L-) (H1045)	1805027701	1
Hepart AG	Isoleucine (L-) (H1045)	15001640/2	1
Hepart AG	L-5-HTP (H1035)	16000441/0	1
Hepart AG	L-5-HTP (H1035)	15001503	1
Hepart AG	L-5-HTP (H1035)	1837999801	1
Hepart AG	L-5-HTP (H1035)	15001503/2	3
Hepart AG	Lavender flower	A130214	4
Hepart AG	Leucine (L-) (H1046)	18001305/0	1
Hepart AG	Leucine (L-) (H1046)	15000283/1	2
Hepart AG	Leucine (L-) (H1046)	1809994101	1
Hepart AG	Lipoic acid (a-) (H10...	14000497/1	1
Hepart AG	Lipoic acid (a-) (H10...	150001526/2	1
Hepart AG	Lipoic acid (a-) (H10...	1313025701	2
Hepart AG	Lipoic acid (a-) (H10...	15001526/2	3
Hepart AG	Lysine (L-) (H1047)	15000372/2	1
Hepart AG	Lysine (L-) (H1047)	16001395/1	3
Hepart AG	Lysine (L-) (H1047)	16001503/0	2
Hepart AG	Lysine (L-) (H1047)	65401	2
Hepart AG	Lysine (L-) (H1047)	16001503/1	1
Hepart AG	Magnesium (H1060)	16001703/0	1
Hepart AG	Magnesium (H1060)	14001895/2	1
Hepart AG	Magnesium (H1060)	16000272/0	2
Hepart AG	Magnesium (H1060)	16001705/0	2
Hepart AG	Magnesium (H1060)	43602	5
Hepart AG	Magnesium (H1060)	43701	4
Hepart AG	Magnesium (H1060)	14001894/0	1
Hepart AG	Magnesium (H1060)	16000269	2
Hepart AG	Magnesium (H1060)	69901	3
Hepart AG	Matricaria flower	B-477, kBA	1
Hepart AG	Mullein flower, cut	12170904	2
Hepart AG	Nettle leaf, cut	P234605kBA	1
Hepart AG	OPC Grape seed (H1066)	16000637	2
Hepart AG	OPC Grape seed (H1066)	1504018401	2
Hepart AG	OPC Grape seed (H1066)	70401	5
Hepart AG	Peppermint leaf, cut	B-185	5
Hepart AG	Phenylalanine (L-) (H...	74401/1	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart AG	Phenylalanine (L-) (H...	17000876/0	1
Hepart AG	Potassium (H1030)	1124944201	2
Hepart AG	Potassium (H1030)	1124028903	1
Hepart AG	Potassium (H1030)	16001645/0	4
Hepart AG	Potassium (H1030)	16001646/0	5
Hepart AG	Proline (-L) (H1051)	25707	2
Hepart AG	Proline (-L) (H1051)	707801	2
Hepart AG	Restharrow root, cut	B130404	3
Hepart AG	Rhodiola (H1064)	15001045/4	1
Hepart AG	Rhodiola (H1064)	70501	1
Hepart AG	Ribwort plantain, cut	20832	5
Hepart AG	Savory, rubbed	12177205	3
Hepart AG	Selenium (H1067)	16001350/0	15
Hepart AG	Selenium (H1067)	933301	1
Hepart AG	Selenium (H1067)	15000157/1	1
Hepart AG	Selenium (H1067)	16001349/0	2
Hepart AG	Silymarin (H1071)	1244029901	1
Hepart AG	Silymarin (H1071)	67701	4
Hepart AG	St. John's wort, cut	B-445,kbA	1
Hepart AG	Stemless carline root	B120726	1
Hepart AG	Taurine (H1084)	16000947/0	1
Hepart AG	Taurine (H1084)	2656905/0	2
Hepart AG	Taurine (H1084)	16000947/2	2
Hepart AG	Thyme, rubbed	A130109	2
Hepart AG	Trace elements JK (H1...	17000870/0	3
Hepart AG	Trace elements JK (H1...	15000528/0	1
Hepart AG	Trace elements JK (H1...	15000525	1
Hepart AG	Trace elements JK (H1...	72701	3
Hepart AG	Trace elements JK (H1...	1026011701	4
Hepart AG	Trace elements SE (H1...	14000720	3
Hepart AG	Trace elements SE (H1...	42301	4
Hepart AG	Trace elements SE (H1...	31301	3
Hepart AG	Trace elements SE (H1...	937201	7
Hepart AG	Trace elements SE (H1...	15000556	1
Hepart AG	Trace elements SE (H1...	41401	1
Hepart AG	Trace elements SE (H1...	120435-14	1
Hepart AG	Tryptophan (L-) (H1053)	27901	7
Hepart AG	Tryptophan (L-) (H1053)	16001404/3	3
Hepart AG	Tryptophan (L-) (H1053)	160001404/3	1
Hepart AG	Tryptophan (L-) (H1053)	1824027901	2
Hepart AG	Tyrosine (L-) (H1054)	74601	2
Hepart AG	Tyrosine (L-) (H1054)	15001380/0	2
Hepart AG	Valine (L-) (H1056)	150014360	1
Hepart AG	Valine (L-) (H1056)	16001380/0	1
Hepart AG	Valine (L-) (H1056)	1807025012	1
Hepart AG	Vitamin B1 (H1101)	16000698/0	1
Hepart AG	Vitamin B1 (H1101)	14001112	1
Hepart AG	Vitamin B1 (H1101)	1622030202	2
Hepart AG	Vitamin B1 (H1101)	53301	3
Hepart AG	Vitamin B1 (H1101)	1621030202	1
Hepart AG	Vitamin B12 (H1102)	1064025017	4
Hepart AG	Vitamin B12 (H1102)	17000276/0	1
Hepart AG	Vitamin B12 (H1102)	1064026401	2
Hepart AG	Vitamin B6 (H1103)	1552025711	2
Hepart AG	Vitamin B6 (H1103)	1554020101	1
Hepart AG	Vitamin C (H1104)	16000197/0	1
Hepart AG	Vitamin C (H1104)	74201	2
Hepart AG	Vitamin C (H1104)	1084025020	4
Hepart AG	Vitamin C (H1104)	44601	2
Hepart AG	Vitamin C (H1104)	44501	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart AG	Vitamin C (H1104)	17000917/0	1
Hepart AG	Vitamin C (H1104)	17000915/0	2
Hepart AG	Vitamin complex (H1100)	1014025301	1
Hepart AG	Vitamin complex (H1100)	16000957/0	1
Hepart AG	Vitamin complex (H1100)	17000880/0	4
Hepart AG	Vitamin complex (H1100)	H110031	1
Hepart AG	Vitamin complex (H1100)	15001277	1
Hepart AG	Vitamin complex (H1100)	49202	6
Hepart AG	Vitamin complex (H1100)	15001275/1	1
Hepart AG	Vitamin complex (H1100)	77201	6
Hepart AG	Vitamin D3 (H1105)	64801	10
Hepart AG	Vitamin D3 (H1105)	15000490/2	1
Hepart AG	Vitamin D3 (H1105)	16000952/1	1
Hepart AG	Vitamin E NAT (H1107)	1094025022	1
Hepart AG	Vitamin E NAT (H1107)	15001557	1
Hepart AG	Vitamin E NAT (H1107)	17000890/0	3
Hepart AG	Vitamin E NAT (H1107)	16001246/0	1
Hepart AG	Vitamin E NAT (H1107)	65101	3
Hepart AG	Vitamin K2-Ca-Mn-Si c...	15000836/6	1
Hepart AG	Walnut leaf, cut	A120123	1
Hepart AG	Willowherb, small flo...	10326110	2
Hepart AG	Zinc (H1110)	16000519/0	10
Hepart AG	Zinc (H1110)	41501	12
Hepart AG	Zinc (H1110)	16000518/1	1
Hepart AG, CH-8280 Kr...	Acacia gum (H1001)	1211622002	2
Hepart AG, CH-8280 Kr...	Antioxidants NAT comp...	52101	4
Hepart AG, CH-8280 Kr...	Antioxidants NAT comp...	1034019901	1
Hepart AG, CH-8280 Kr...	Antioxidants NAT comp...	1034030901	2
Hepart AG, CH-8280 Kr...	Betamethasone valerate	11131914	3
Hepart AG, CH-8280 Kr...	Bioflavonoid complex ...	56001	3
Hepart AG, CH-8280 Kr...	Biotin (H1007)	1592029603	4
Hepart AG, CH-8280 Kr...	Calcium (H1008)	1104008901	10
Hepart AG, CH-8280 Kr...	Carnitine (L-) (H1038)	26020	2
Hepart AG, CH-8280 Kr...	Carnitine (L-) (H1038)	1184965901	1
Hepart AG, CH-8280 Kr...	Catechin extract (H10...	62901	2
Hepart AG, CH-8280 Kr...	Chondroitin sulfate ...	1820031801	1
Hepart AG, CH-8280 Kr...	Glutathione (H1020)	1251025703	1
Hepart AG, CH-8280 Kr...	Glutathione (H1020)	1251019703	2
Hepart AG, CH-8280 Kr...	Iron (H1015)	1114998701	5
Hepart AG, CH-8280 Kr...	Isoleucine (L-) (H1045)	66301	3
Hepart AG, CH-8280 Kr...	Isoleucine (L-) (H1045)	1804968701	1
Hepart AG, CH-8280 Kr...	Lysine (L-) (H1047)	1800955002	1
Hepart AG, CH-8280 Kr...	Magnesium (H1060)	1134973501	9
Hepart AG, CH-8280 Kr...	Magnesium (H1060)	1134043601	1
Hepart AG, CH-8280 Kr...	OPC Grape seed (H1066)	130208	1
Hepart AG, CH-8280 Kr...	Silymarin (H1071)	1241046101	1
Hepart AG, CH-8280 Kr...	Silymarin (H1071)	1241950401	1
Hepart AG, CH-8280 Kr...	Taurine (H1084)	25710	3
Hepart AG, CH-8280 Kr...	Vitamin B6 (H1103)	25711	5
Hepart AG, CH-8280 Kr...	Vitamin C (H1104)	1084025019	1
Hepart AG, CH-8280 Kr...	Vitamin complex (H1100)	49301	14
Hepart AG, CH-8280 Kr...	Vitamin complex (H1100)	1014025101	2
Hepart AG, CH-8280 Kr...	Vitamin D3 (H1105)	1336008101	2
Hepart AG, CH-8280 Kr...	Vitamin E NAT (H1107)	65001	4
Hepart AG, Unisan GmbH	Vitamin B6 (H1103)	15000279/0	2
Hepart/Unisan	Antioxidants NAT comp...	16000376/0	1
Hepart/Unisan	Antioxidants NAT comp...	16000967	1
Hepart/Unisan	Antioxidants NAT comp...	15000877/1	1
Hepart/Unisan	Calcium (H1008)	14001938/1	1
Hepart/Unisan	Folic acid (H1016)	15000125/0	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Hepart/Unisan	Iron (H1015)	15001258/1	2
Hepart/Unisan	Magnesium (H1060)	16000267/0	1
Hepart/Unisan	Taurine (H1084)	14001767/3	5
Hepart/Unisan	Threonine (L-) (H1052)	14001718/3	1
Hepart/Unisan	Zinc (H1110)	14001803/2	2
Herba	Cetostearyl alcohol	KNC1133J18	1
Herbal creations/Caes...	Methoxsalen	15244208	1
Herbal creations/Caes...	Methoxsalen	15244212	1
Herbal Creations/Caes...	Methoxsalen	15244203	1
Herbal Creations/Caes...	Methoxsalen	15244205	3
Herbasin	Huang Qin	g225h120442	1
Herbasin	Lai Fu Zi	g209h1220321	1
Herbasinica	(Shi) Chang Pu	g006h0544321	1
Herbasinica	Agrimony, cut	11374007	1
Herbasinica	Lavender flower	1488A130830	1
Herstellung	Captopril	12339607	1
Hetterich/Noewda	Neostigmine bromide	17131901	2
HiperScan	Peppermint leaf, cut		1
Hirsch	Alfatradiol	12L19-N08	9
Hirsch	Betamethasone valerate	13G19-No1	1
Hirsch	Tannin	1811209	1
Hirsch	Urea	810811	1
Hirundo Products Vert...	Arginine hydrochloride	6412784-152	2
Hof-Apotheke,St.Afra	Comfrey root	125041	1
Holdermann	Betamethasone valerate	120801402	5
Holdermann	Mometasone furoate, m...	14C07-B05-29126	1
Holdermann	Salicylic acid	14E26-B04294624	1
Hubei Hongyuan	Metronidazole	140700301	3
Hubertus Apotheke	Equisetum stem, cut	3831C140403	1
Ichthyol	Wax, yellow	13031239	12
Ichthyol Gesellschaft...	Fluocinolonacetoneid, ...	14L19-B04-321003	1
Ichthyol Gesellschaft...	Biotin (H1007)	1591983601	2
Ichthyol/Noweda	Minoxidil	17117803	1
Ichthyol/Phoenix	Fenugreek, powdered	14377205	1
Ichthyol/Phoenix	Miconazole nitrate	16/098/M	1
Ichthyol/Phoenix	Sorbitan stearate	14372803	2
Ichthyol/Phönix	Carbomers 35000 / 50000	14149905	3
Ichthyol/Phönix	Magnesium hydrogen ci...	14/12-B02-298854	3
Ichthyol/Phönix	Sodium dodecyl sulfate	15344510	6
Ichthyol-Gesellschaft	Betamethasone valerate	13k29-B32-288247	1
Ichthyol-Gesellschaft	Erythromycin	161020E/1607038-01	1
ICHTHYOL-GESELLSCHAFT	Sodium monohydrogenph...	Caelo	1
Ichthyol Gesellschaft...	Citric acid	16038618	1
Ichthyol/Gehe	Triamcinolone acetoneid	15430102	5
Ichthyol/Phönix	Alchemilla, cut	13141107	3
Ichthyol/Phönix	Zinc orotate dihydrate	12B02-N02	1
ICROM	Minoxidil	6/12	3
icrom spa	Minoxidil	511v	1
Imhof	Menthol	15-5401	1
inalco	D-galactose	212619PF5	1
inalco	D-galactose	42619PF4	1
inalco	D-galactose	32619PF8	1
Industrias GmbH / Ph	Sulfacetamide	64518	1
InfectoPharm	Propranolol hydrochlo...	J0512011	1
inresa	Amifampridine	58710001	2
inresa	Amifampridine	58010001	1
Inresa	Amifampridine	33-16-I-77P3018/1307	2
Inresa	Amifampridine	33-15-IV-56P2749B/22	3
Inresa	Amifampridine	33-16-I-77P3018	2
Inresa	Amifampridine	580710001	11

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Inresa	Amifampridine	631350001	4
Inresa	Amifampridine	63150001	1
Inresa	Amifampridine	508670001	1
Inresa	Amifampridine	33-15-IV-56P2749B	2
Inresa	Amifampridine	33-15-IV-56P2749	2
Inresa	Fampridine	220825	1
Inresa	Melatonin	110005	1
Inresa	Metoprolol tartrate	504005RIIC/29082016	1
Inresa	Metoprolol tartrate	4050005RIIB	1
Inresa 22.09.15	Propranolol hydrochloride	3043PRRII	1
Inresa Arznei	Pyridoxal-5-phosphate	13100411	1
Intern Drogen	Sodium dodecyl sulfate	216991080	1
INTR/ Sana	Soy lecithin (instant...)	16223302	3
IPCA/Noweda	Furosemide	11D29-N08	1
Jenapharm/Phönix	Erythromycin	1607036-01	1
Jenne	Aniseed	2700A131015	2
Jenne	Carbomers 35000 / 50000	18049301	1
Jenne	Citric acid	15063002	1
Jenne	Clobetasol propionate	17325601	8
Jenne	Fructose	14A16-B07	1
Jenne	Glucose	12141705	2
Jenne	Paracetamol	15245901	3
Jenne	Phenylephrine hydrochloride	14H13-B02	2
Jenne	Saccharin sodium	15087707	1
Jenne	Salicylic acid	12345	1
Jenne	Sodium monohydrogenphosphate	12109003	3
Jenne	Tetracaine hydrochloride	13D11-No7	1
Jenne	Trometamol	17320811	21
Jenne	Vanillin	18073501	4
JENNE	Tragacanth, powdered	16110801	1
Jungbunzlauer Ladenbau...	Citric acid	1156158	2
k	Matricaria flower	1290A140513	1
Kaufland	Sucrose	L24494	2
Kaufland	Sucrose	L24475	1
KBA Ware	Equisetum stem, cut	21879	6
kehr	Sodium edetate	17D25-B08-339881	4
Kehr	alpha-tocopheryl acetate	1705003-03	3
Kehr	Ambroxol hydrochloride	18013201	6
Kehr	Benzocaine	15C19-B04-325359	3
Kehr	Benzocaine	16D13-B01-330168	2
Kehr	Cetyl alcohol	17047303	1
Kehr	Diclofenac sodium	1711026-03	6
Kehr	Gentamicin sulfate	110036-02	1
Kehr	Hydroxyethylcellulose	211112PML01	1
Kehr	Lidocaine	17G03-B08-338284	1
Kehr	Lidocaine	16F21-B06-328421	10
Kehr	Lidocaine	16F21-B06-332251	1
Kehr	Nicotinic acid	1801021-01	2
Kehr	Nystatin	1612001-02	1
Kehr	Oxetacaine	17249805	1
Kehr	Sodium benzoate	1701045-01	2
Kehr	Stearic acid 50	16346902	1
Kehr	Triamcinolone acetonide	1057902	1
Kehr	Urea	1802007	1
keine Angabe	Chlorhexidine diacetate	11204402	2
Kenk / Phoenix	Artichoke leaf, cut	1981B130111	3
Kenk / Phoenix	Daisy flower	1210A130307	4
Kenk / Phoenix	European golden rod, ...	1781B110826	1
Kenk / Phoenix	Fumitory, cut	3881B111013	2
Kenk / Phoenix	Hawthorn leaf with fl. ...	1971A130307	6

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Kenk / Phoenix	Lapacho wood, cut	5161A110204	1
Kenk / Phoenix	Raspberry leaf, cut	22781C120508	1
Kenk / Phoenix	Sage leaf, cut	2301A130904	2
Kenk / Phoenix	Turmeric, javanese, cut	6731B130117	1
Klenk	Abrotani herb, cut	21081301	1
Klenk	Abrotani herb, cut	26011403	1
Klenk	Acorus root, natural,...	B130807	1
Klenk	Acorus root, peeled	6701B150407	1
Klenk	Acorus root, peeled	6701B130870	1
Klenk	Acorus root, peeled	6701b130103	1
Klenk	Acorus root, peeled	6701B120817	6
Klenk	Acorus root, peeled	6701B120301	3
Klenk	Acorus root, peeled	6701b120301	1
Klenk	Acorus root, peeled	6701B120427	7
Klenk	Acorus root, peeled	6711B110606	1
Klenk	Acorus root, peeled	106212173	1
Klenk	Acorus root, peeled	4797362	1
Klenk	Agnus castus fruit	A120803	1
Klenk	Agnus castus fruit	A110614	3
Klenk	Agnus castus fruit	12274007	3
Klenk	Agnus castus fruit	7110A120803	1
Klenk	Agnus castus fruit	7110A110614	2
Klenk	Agnus castus fruit	A130613	1
Klenk	Agnus castus fruit	7110N150121	1
Klenk	Agnus castus fruit	7110A130613	1
Klenk	Agnus castus fruit	7110N090701	1
Klenk	Agrimony, cut	3481C120508	1
Klenk	Alchemilla, cut	3501B090417	1
Klenk	Alchemilla, cut	3501A110301	3
Klenk	Alchemilla, cut	3501C121105	4
Klenk	Alchemilla, cut	3501B100428	1
Klenk	Almonds, bitter, whole	40A111004	1
Klenk	Almonds, bitter, whole	40A120924	4
Klenk	Almonds, bitter, whole	40A121010	11
Klenk	Almonds, bitter, whole	40A130909	9
Klenk	Almonds, bitter, whole	40A140908	19
Klenk	Almonds, bitter, whole	40A0130909	1
Klenk	Almonds, bitter, whole	40A111117	7
Klenk	Almonds, bitter, whole	40A150904	12
Klenk	Almonds, bitter, whole	40A161028	12
Klenk	Almonds, bitter, whole	40A150905	1
Klenk	Almonds, bitter, whole	40A0161005	1
Klenk	Almonds, bitter, whole	40A161206	8
Klenk	Almonds, bitter, whole	40A-150904	2
Klenk	Angelica root	6081N120110	14
Klenk	Angelica root	6081b130109	2
Klenk	Angelica root	6081B130410	4
Klenk	Angelica root	6082b130410	1
Klenk	Angelica root	6081B170601	1
Klenk	Angelica root	6091B130109	1
Klenk	Angelica root	B130109	2
Klenk	Angelica root	6081B140918	1
Klenk	Angelica root	6081n120110	1
Klenk	Angelica root	6081B130726	1
Klenk	Angelica root	6081B150701	3
Klenk	Angelica root	6081B130109	1
Klenk	Angelica root		1
Klenk	Angelica root	6085B120614	1
Klenk	Aniseed	212M-01530	1
Klenk	Aniseed	814Q01530	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Aniseed	2700A-110701	1
Klenk	Aniseed	2700A160511	5
Klenk	Aniseed	2700A150601	5
Klenk	Aniseed	2700A110701	7
Klenk	Aniseed	2700A120320	8
Klenk	Aniseed	2901311530	1
Klenk	Aniseed	2103311530	1
Klenk	Aniseed	3113Q01530	1
Klenk	Aniseed	2700a120320	1
Klenk	Aniseed	2700A121107	1
Klenk	Aniseed	2700A130527	13
Klenk	Aniseed	2700A12107	1
Klenk	Aniseed	2700A150901	1
Klenk	Aniseed	2700A140129	6
Klenk	Aniseed, pounded	2706B100120	1
Klenk	Aniseed, powdered	2704B110824	1
Klenk	Apple, cut	A120302	3
Klenk	Arnica flower	1181B130605	1
Klenk	Artichoke leaf, cut	B121025	2
Klenk	Artichoke leaf, cut	1981B121025	2
Klenk	Artichoke leaf, cut	1981B120413	1
Klenk	Ash leaf, cut	5111A-01491	1
Klenk	Asparagus root, cut	6151B100332	1
Klenk	Asplenium scolopendri...	4451B091116	1
Klenk	Asplenium scolopendri...	4451B120605	1
Klenk	Barberry root bark, cut	551B120320	1
Klenk	Barberry root bark, cut	551B110128	4
Klenk	Basil, rubbed	A120502	1
Klenk	Basil, rubbed	3632a120502	1
Klenk	Basil, rubbed	3632A120502	4
Klenk	Basil, rubbed	3632A110201	3
Klenk	Basil, rubbed	A110201	1
Klenk	Bean shells, cut	671A100816	2
Klenk	Bean shells, cut	671A110301	1
Klenk	Bean shells, cut	171A110301	1
Klenk	Bean shells, cut	671C120321	16
Klenk	Bean shells, cut	671A-110420	1
Klenk	Bean shells, cut	671A110420	9
Klenk	Bean shells, cut	C121105	2
Klenk	Bean shells, cut	C120321	2
Klenk	Bean shells, cut	671C130620	6
Klenk	Bean shells, cut	C130620	1
Klenk	Bean shells, cut	671C121105	8
Klenk	Bearberry leaf, cut	2391B07030901	1
Klenk	Bearberry leaf, cut	2391B130328	9
Klenk	Bearberry leaf, cut	2391B120719	13
Klenk	Bearberry leaf, cut	2391B120314	6
Klenk	Bearberry leaf, cut	2391b120314	2
Klenk	Bilberry fruit, dried	910N130131	1
Klenk	Bilberry fruit, dried	2910N121022	7
Klenk	Bilberry fruit, dried	N121022	1
Klenk	Bilberry fruit, dried	2910N130131	11
Klenk	Bilberry fruit, dried	2910n130131	1
Klenk	Bilberry fruit, dried	2910A090717	1
Klenk	Bilberry fruit, dried	2010N110304	1
Klenk	Bilberry fruit, dried	N130131	2
Klenk	Bilberry fruit, dried	2910n110304	1
Klenk	Bilberry fruit, dried	2910N110304	5
Klenk	Bilberry leaf, cut	N110516	1
Klenk	Bilberry leaf, cut	4042013	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Bilberry leaf, cut	2191N110516	5
Klenk	Birch leaf, cut	A121130	2
Klenk	Birch leaf, cut	1901C120802	19
Klenk	Birch leaf, cut	F-040413-1	1
Klenk	Birch leaf, cut	1901A130318	2
Klenk	Birch leaf, cut	C120802	1
Klenk	Birch leaf, cut	21061303	1
Klenk	Birch leaf, cut	1901A140214	1
Klenk	Birch leaf, cut	4113A-01481	1
Klenk	Birch leaf, cut	1901c120802	1
Klenk	Birch leaf, cut	4021409	1
Klenk	Birch leaf, cut	1901A05051801	1
Klenk	Birch leaf, cut	20111201	1
Klenk	Birch leaf, cut	4213I01481	1
Klenk	Birch leaf, cut	28101306	1
Klenk	Birch leaf, cut	1901c120508	1
Klenk	Birch leaf, cut	1901A-110726	3
Klenk	Birch leaf, cut	27121301	1
Klenk	Birch leaf, cut	1901A110726	24
Klenk	Birch leaf, cut	1901A130218	4
Klenk	Birch leaf, cut	1901A121130	24
Klenk	Birch leaf, cut	1901C120508	18
Klenk	Bismuth gallath, basic	12E16-N03	3
Klenk	Bitter fennel	2866B100211	1
Klenk	Bitter fennel	2880N120418	4
Klenk	Bitter fennel	2860A130226	6
Klenk	Bitter fennel	2860a120412	3
Klenk	Bitter fennel	2860A150320	12
Klenk	Bitter fennel	3812A-01545	1
Klenk	Bitter fennel	2860A150617	4
Klenk	Bitter fennel	28060A140612	1
Klenk	Bitter fennel	2860A140613	2
Klenk	Bitter fennel	2880A100503	2
Klenk	Bitter fennel	2860A140612	12
Klenk	Bitter fennel	2860A140224	1
Klenk	Bitter fennel	2860A140224	1
Klenk	Bitter fennel	2860, A150617	1
Klenk	Bitter fennel	2866B110318	1
Klenk	Bitter fennel	2866B130206	1
Klenk	Bitter fennel	2860A110704	4
Klenk	Bitter fennel	2860A130513	2
Klenk	Bitter orange flower ...	3613m-01445	1
Klenk	Bitter orange flower ...	1201B110809	1
Klenk	Bitter orange flower ...	1201B120822	1
Klenk	Bitter orange peel, cut	511B120830	4
Klenk	Bitter orange peel, cut	3613M-01531	1
Klenk	Bitter orange peel, cut	B120830	1
Klenk	Bitter orange peel, cut	2711B111207	1
Klenk	Bitter orange peel, cut	511C130226	1
Klenk	Bitter orange peel, cut	511B110321	2
Klenk	Bitter orange peel, cut	511A110411	1
Klenk	Bitter orange peel, cut	C130226	3
Klenk	Black cohosh, cut	6721B141202	1
Klenk	Black cohosh, cut	3794377	1
Klenk	Black cohosh, cut	6721B130214	10
Klenk	Black cohosh, cut	6721B081208	1
Klenk	Black cohosh, cut	6721B1130528	1
Klenk	Black cohosh, cut	6721b130528	1
Klenk	Black cohosh, cut	6721B130528	3
Klenk	Black cohosh, cut	B130528	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Black cumin seed	A120703	1
Klenk	Black cumin seed	7330A120703	1
Klenk	Black mustard seed	7430A120813	2
Klenk	Black mustard seed, p...	7434A120712	1
Klenk	Black mustard seed, p...	7434A130422	3
Klenk	Black mustard seed, p...	7434A140923	2
Klenk	Black mustard seed, p...	A130422	1
Klenk	Black mustard seed, p...	7434A150608	4
Klenk	Black mustard seed, p...	7434a150608	1
Klenk	Black pepper	A130415	2
Klenk	Blackberry leaf, cut	2271B111216	1
Klenk	Blackcurrant leaf, cut	2251a110204	1
Klenk	Blackcurrants	2990A110620	1
Klenk	Blackcurrants	A121106	1
Klenk	Blackthorn flower	2612E-01468	1
Klenk	Blackthorn flower, ru...	1150n101109	1
Klenk	Blackthorn flower, ru...	N101109	2
Klenk	Blessed thistle, cut	A110517	1
Klenk	Blue bottle flower wi...		1
Klenk	Boldo leaf, cut	1911A100714	7
Klenk	Boldo leaf, cut	1911B090709	1
Klenk	Buckwheat herb, cut	4581B120509	3
Klenk	Burdock root, cut	A120206	1
Klenk	Burdock root, cut	6161A120208	4
Klenk	Burdock root, cut	6161A-120103	1
Klenk	Burdock root, cut	6161A080709	1
Klenk	Burdock root, cut	6161A120103	1
Klenk	Burdock root, cut	6161A101206	1
Klenk	Burdock root, cut	6161A-101206	1
Klenk	Burdock root, cut	A130214	1
Klenk	Calendula flower	1241B130517	1
Klenk	Calendula flower	1241B120312	1
Klenk	Calendula flower	C120321	1
Klenk	Calendula flower	1250A110511	1
Klenk	Calendula flower, cut	1241B130927	1
Klenk	Calendula flower, cut	1241B160530	1
Klenk	Calendula flower, cut	1241B110118	2
Klenk	Calendula flower, cut	1241B120419	2
Klenk	Calendula flower, cut	1241B130517	5
Klenk	Calendula flower, cut	112211448	1
Klenk	Calendula flower, cut	1240A110418	1
Klenk	Calendula flower, cut	1241N151209	1
Klenk	Calendula flower, cut	204311448	1
Klenk	Calendula flower, cut	1241B130321	5
Klenk	Calendula flower, cut	1241B140211	1
Klenk	Calendula flower, cut	B130517	2
Klenk	Calendula flower, cut	1241C120321	4
Klenk	Calendula flower, cut	3013E1448	1
Klenk	Calendula flower, cut	1241B150902	1
Klenk	Camphor	444A160817	6
Klenk	Camphor	444A140526	4
Klenk	Camphor	444A100826	1
Klenk	Camphor	444A150721	1
Klenk	Camphor	444A150721	1
Klenk	Caraway fruit	5213A1537	1
Klenk	Caraway fruit	2770A131008	2
Klenk	Caraway fruit	2770A130910	1
Klenk	Caraway fruit	2770A160304	1
Klenk	Caraway fruit	412091-1-P	1
Klenk	Caraway fruit	2770A110601	9

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Caraway fruit	2770A140708	4
Klenk	Caraway fruit	2770A110818	8
Klenk	Caraway fruit	2770A130521	12
Klenk	Caraway fruit	2770A120403	14
Klenk	Caraway fruit	1810A-01537	1
Klenk	Caraway fruit	A120905	3
Klenk	Caraway fruit	2770a120905	2
Klenk	Caraway fruit	2270A130521	1
Klenk	Cardamom fruit	A110810	3
Klenk	Cardamom fruit	2750A110810	2
Klenk	Cardamom fruit	A130807	1
Klenk	Centaury, cut	3711B-111216	1
Klenk	Centaury, cut	B120628	1
Klenk	Centaury, cut	3711B111216	1
Klenk	Ceylon cinnamon, fine...	615B130325	3
Klenk	Ceylon cinnamon, fine...	611B100628	1
Klenk	Ceylon cinnamon, fine...	614B121024	4
Klenk	Ceylon cinnamon, fine...	614B121010	2
Klenk	Ceylon cinnamon, fine...	614B120411	2
Klenk	Ceylon cinnamon, fine...	614b121024	1
Klenk	Ceylon cinnamon, fine...	614a120411	1
Klenk	Chelidonium majus her...	3731A120611	4
Klenk	Chelidonium majus her...	3731K120904	4
Klenk	Chelidonium majus her...	3731A110301	6
Klenk	Chelidonium majus her...	1406111632	1
Klenk	Chelidonium majus her...	3731B100205	1
Klenk	Chelidonium majus her...	4711E-01632	1
Klenk	Chelidonium majus her...	3731A130904	2
Klenk	Chickweed, cut	4551b110222	1
Klenk	Chickweed, cut	4551B090528	1
Klenk	Chicory herb, cut	3751A120307	7
Klenk	Chicory herb, cut	3751A120123	2
Klenk	Chicory herb, cut	3751A-100118	1
Klenk	Chicory herb, cut	3751c120614	1
Klenk	Chicory herb, cut	3751a120307	1
Klenk	Chicory herb, cut	4510I-01634	1
Klenk	Chicory herb, cut	A120123	2
Klenk	Chicory herb, cut	3751A080902	1
Klenk	Chicory root, cut	66221A101202	1
Klenk	Chicory root, cut	6221A130304	1
Klenk	Chicory root, cut	6221A090311	1
Klenk	Chicory root, cut	6221N130201	1
Klenk	Chicory root, cut	A130304	1
Klenk	Cinchona bark, cut	601B130415	1
Klenk	Cinchona bark, cut	601B130415	1
Klenk	Cinnamon rods (Ceylon...	610A101018	1
Klenk	Cinnamon rods (Ceylon...	B121016	1
Klenk	Cinnamon rods (Ceylon...	A101018	1
Klenk	Cinnamon, cut	610A130927	1
Klenk	Cinnamon, cut	610A101018	1
Klenk	Clove	A120328	1
Klenk	Clove	1260A100617	2
Klenk	Clove	1260A150304	2
Klenk	Clove	1260A130903	3
Klenk	Clove	1260A120328	5
Klenk	Clove	1260A101223	1
Klenk	Clove	1260A110608	5
Klenk	Coltsfoot leaf, cut	2021C140730	1
Klenk	Coltsfoot leaf, cut	2021B150407	1
Klenk	Coltsfoot leaf, cut	23121302	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Coltsfoot leaf, cut	5122013a	1
Klenk	Coltsfoot leaf, cut	2021B140131	1
Klenk	Coltsfoot leaf, cut	2021B111005	1
Klenk	Coltsfoot leaf, cut	B130227	1
Klenk	Coltsfoot leaf, cut	2021B130319	4
Klenk	Coltsfoot leaf, cut	20021301	1
Klenk	Coltsfoot leaf, cut	14111302	1
Klenk	Coltsfoot leaf, cut	5011301	1
Klenk	Coltsfoot leaf, cut	11101202	1
Klenk	Coltsfoot leaf, cut	2021B121203	5
Klenk	Coltsfoot leaf, cut	26091301	1
Klenk	Coltsfoot leaf, cut	2021C130326	8
Klenk	Coltsfoot leaf, cut	2021B120124	4
Klenk	Coltsfoot leaf, cut	2020B120124	1
Klenk	Coltsfoot leaf, cut	8061207	1
Klenk	Coltsfoot leaf, cut	2021b121204	1
Klenk	Coltsfoot leaf, cut	2021B120601	5
Klenk	Coltsfoot leaf, cut	KS100061	1
Klenk	Coltsfoot leaf, cut	2021B130227	12
Klenk	Coltsfoot leaf, cut	B120124	1
Klenk	Comfrey herb, cut	4561B100518	2
Klenk	Comfrey root	6241B120531	3
Klenk	Comfrey root	6241b130724	1
Klenk	Comfrey root	6241B131009	1
Klenk	Comfrey root	6241B130724	7
Klenk	Comfrey root	6241b130619	1
Klenk	Comfrey root	6244B130624	1
Klenk	Comfrey root	6241b120807	1
Klenk	Comfrey root	B130619	1
Klenk	Comfrey root	6241B130619	8
Klenk	Comfrey root	6241B120807	4
Klenk	Comfrey root	6241B130312	3
Klenk	Comfrey root	6241B100715	1
Klenk	Comfrey root	6241B06110202	1
Klenk	Comfrey root	6241B121214	2
Klenk	Common wormwood herb,...	2591A130218	1
Klenk	Common wormwood herb,...	3591A130218	10
Klenk	Common wormwood herb,...	3591A120123	8
Klenk	Common wormwood herb,...	3591a130218	2
Klenk	Common wormwood herb,...	A120123	2
Klenk	Common wormwood herb,...	3591A11119	1
Klenk	Common wormwood herb,...	3591a120123	1
Klenk	Coriander	2790A161111	5
Klenk	Coriander	2709A130129	1
Klenk	Coriander	2790A130129	2
Klenk	Coriander	2790a121023	1
Klenk	Coriander	2790A121023	7
Klenk	Coriander	A121023	3
Klenk	Coriander	2790A151117	1
Klenk	Corn starch	64A130304	10
Klenk	Corn starch	64A170123	6
Klenk	Corn starch	64A080218	4
Klenk	Corn starch	4263535	1
Klenk	Cottonwood bud	3260A120828	1
Klenk	Cowberry leaf, cut	2431B101006	2
Klenk	Daisy flower	1210a071301	1
Klenk	Daisy flower	1210A110223	1
Klenk	Daisy flower	A130705	1
Klenk	Daisy flower	1210A120215	3
Klenk	Damiana leaf, cut	A130206	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Damiana leaf, cut	1990A130206	1
Klenk	Damiana leaf, cut	12064101	1
Klenk	Dandelion herb with r...	6611B1104104	1
Klenk	Dandelion herb with r...	6611B-131206	1
Klenk	Dandelion herb with r...	5122013b	1
Klenk	Dandelion herb with r...	6611B090306	1
Klenk	Dandelion herb with r...	6611B-120620	1
Klenk	Dandelion herb with r...	6611B130314	2
Klenk	Dandelion herb with r...	6611B110929	2
Klenk	Dandelion herb with r...	6611A110104	1
Klenk	Dandelion herb with r...	6611B120620	19
Klenk	Dandelion herb with r...	6611B110131	1
Klenk	Dandelion herb with r...	6611B110104	5
Klenk	Dandelion herb with r...	6611B160613	1
Klenk	Dandelion root, cut	6621N130813	6
Klenk	Dandelion root, cut	N130813	1
Klenk	Dandelion root, cut	A120606	1
Klenk	Dandelion root, cut	6621A120606	15
Klenk	Dandelion root, cut	6621A110420	1
Klenk	Dandelion root, cut	6621A141027	5
Klenk	Dandelion root, cut	6621A160111	1
Klenk	Dandelion root, cut	6101202	1
Klenk	Dandelion root, cut	6621A120208	6
Klenk	Dandelion root, cut	6621A140114	3
Klenk	Dandelion root, cut	6621A140114	3
Klenk	Devil's claw root	6331B130129	5
Klenk	Devil's claw root	6331B120508	3
Klenk	Devil's claw root	10111204	1
Klenk	Dill fruit	2680A120312	7
Klenk	Dill fruit	2680a120312	1
Klenk	Dill fruit	2680A121211	3
Klenk	Dill fruit	2680A-120312	1
Klenk	Dill fruit	A120312	1
Klenk	Dill herb, cut	3531B110726	1
Klenk	Elder berries	3010N120921	3
Klenk	Elder flower, rubbed	1571A-110223	1
Klenk	Elder flower, rubbed	16061201	1
Klenk	Elder flower, rubbed	26061201	1
Klenk	Elder flower, rubbed	1571A130627	2
Klenk	Elder flower, rubbed	21091201	1
Klenk	Elder flower, rubbed	1571A110125	2
Klenk	Elder flower, rubbed	1571A120215	4
Klenk	Elder flower, rubbed	13021301	1
Klenk	Elder flower, rubbed	1571a120215	1
Klenk	Elder flower, rubbed	1571A121121	20
Klenk	Elder flower, rubbed	1571A1202015	1
Klenk	Elder flower, rubbed	A121121	1
Klenk	Elder flower, rubbed	A120215	1
Klenk	Elder flower, rubbed	1571A130930	3
Klenk	Elder flower, rubbed	1571A110223	16
Klenk	Elder flower, rubbed	1571A141126	1
Klenk	Elder flower, rubbed	1571a1211	1
Klenk	Elder root, cut	6281B110427	4
Klenk	Elecampane rhizome, cut	6311141B130423	1
Klenk	Elecampane rhizome, cut	6341B120611	2
Klenk	Elecampane rhizome, cut	6341B130423	3
Klenk	Elecampane rhizome, cut	6341B100426	2
Klenk	Elecampane rhizome, cut	B130423	3
Klenk	Eleutherococcus, cut	7011A110714	1
Klenk	Eleutherococcus, cut	7011A101206	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Eleutherococcus, cut	7011A121204	1
Klenk	Elm bark, cut	841a120413	1
Klenk	Elm bark, cut	841A1101120	1
Klenk	Equisetum stem, cut	3831A120903	8
Klenk	Equisetum stem, cut	3831A120018	1
Klenk	Equisetum stem, cut	3831A111124	4
Klenk	Equisetum stem, cut	3831A1230726	1
Klenk	Equisetum stem, cut	3831c131104	1
Klenk	Equisetum stem, cut	4713U01642	1
Klenk	Equisetum stem, cut	3831a120611	1
Klenk	Equisetum stem, cut	3831A110601	1
Klenk	Equisetum stem, cut	2213U01642	1
Klenk	Equisetum stem, cut	3831A120726	18
Klenk	Equisetum stem, cut	3831A121018	6
Klenk	Equisetum stem, cut	4413Q1642	1
Klenk	Equisetum stem, cut	3831A110922	1
Klenk	Equisetum stem, cut	3831A120611	10
Klenk	Equisetum stem, cut	3831A120529	8
Klenk	Equisetum stem, cut	3831C131104	1
Klenk	Equisetum stem, cut	3831A111124	1
Klenk	Erythromycin	EB068I14	1
Klenk	Estradiol benzoate	8(18144402)	1
Klenk	Eucalyptus leaf, cut	2011B130306-35294110	1
Klenk	Eucalyptus leaf, cut	2011B120419	6
Klenk	Eucalyptus leaf, cut	2011B130306	4
Klenk	Eucalyptus leaf, cut	B120419	1
Klenk	Eucalyptus leaf, cut	2011b130306	1
Klenk	Eucalyptus leaf, cut	B130306	2
Klenk	Euphrasy herb, cut	B101206	1
Klenk	Euphrasy herb, cut	3871B130430	1
Klenk	Euphrasy herb, cut	15029806	1
Klenk	Euphrasy herb, cut	B130430	1
Klenk	European golden rod, ...	1781B110706	1
Klenk	Fenugreek	A130603	1
Klenk	Fenugreek	8061301	1
Klenk	Fenugreek	A120118	2
Klenk	Fenugreek	7260A110530	2
Klenk	Fenugreek	7260A120118	2
Klenk	Fenugreek	7260A130603	3
Klenk	Fenugreek	9011301	1
Klenk	Fenugreek, powdered	A110201	1
Klenk	Fenugreek, powdered	12082014A	1
Klenk	Fenugreek, powdered	7264A110201	9
Klenk	Flax weed herb, cut	4011B131115	3
Klenk	Flax weed herb, cut	4011B121207	2
Klenk	Frangula bark, cut	661B130306	1
Klenk	Frangula bark, cut	661B-110412	1
Klenk	Frangula bark, cut	661B-090903	1
Klenk	Frangula bark, cut	6621B110412	1
Klenk	Frangula bark, cut	661B110412	4
Klenk	Fraxinella root, cut	6271B140708	1
Klenk	Fumitory, cut	B120410	1
Klenk	Fumitory, cut	3881B100127	1
Klenk	Galeopsis herb, cut	3911A130110	3
Klenk	Galeopsis herb, cut	3911A090219	2
Klenk	Galeopsis herb, cut	3911A080409	1
Klenk	Galeopsis herb, cut	3911a130320	1
Klenk	Gentian root, cut	1911E-02142	1
Klenk	Gentian root, cut	6311B130205	2
Klenk	Gentian root, cut	6311a120723	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Gentian root, cut	6311B140704	1
Klenk	Gentian root, cut	6311A120723	3
Klenk	Ginger, peeled	6951B100302	1
Klenk	Ginger, peeled	6951B110119	1
Klenk	Ginger, peeled	6951B110822	1
Klenk	Ginger, peeled	6951B130110	1
Klenk	Ginger, peeled, powde...	6954B121030	2
Klenk	Ginger, peeled, powde...	B121030	2
Klenk	Ginger, peeled, powde...	6954B130418	2
Klenk	Ginger, peeled, powde...	6954B110801	1
Klenk	Ginger, peeled, powde...	6954B170127	1
Klenk	Ginger, peeled, powde...	6954B16029	1
Klenk	Ginger, peeled, powde...	6954B160429	5
Klenk	Ginger, peeled, powde...	KS200260	1
Klenk	Ginger, peeled, powde...	6954B160929	2
Klenk	Goat's rue, cut	3901B140801	1
Klenk	Goat's rue, cut	3901C110905	1
Klenk	Ground ivy herb, cut	3991B110304	1
Klenk	Ground ivy herb, cut	3991B130502	2
Klenk	Hamamelis leaf, cut	2051B110214	1
Klenk	Hamamelis leaf, cut	B110214	1
Klenk	Hawthorn flower	1351B111207	4
Klenk	Hawthorn flower	B1112107	1
Klenk	Hawthorn fruit	2806B121203	2
Klenk	Hawthorn fruit	2806B110818	2
Klenk	Hawthorn leaf with fl...	1971a120215	1
Klenk	Hawthorn leaf with fl...	1971A120727	13
Klenk	Hawthorn leaf with fl...	1971A110620	5
Klenk	Hawthorn leaf with fl...	1972A120727	1
Klenk	Hawthorn leaf with fl...	1971A120215	3
Klenk	Hawthorn leaf with fl...	A130307	1
Klenk	Hawthorn leaf with fl...	1971A130307-35014110	1
Klenk	Hawthorn leaf with fl...	1971A130307-3447	1
Klenk	Hawthorn leaf with fl...	1971a120813	2
Klenk	Hawthorn leaf with fl...	1971A120813	20
Klenk	Hawthorn leaf with fl...	1971A120813-34103101	1
Klenk	Hawthorn leaf with fl...	28101307	1
Klenk	Hawthorn leaf with fl...	14101301	1
Klenk	Hawthorn leaf with fl...	20121302	1
Klenk	Hay flower	10192103	1
Klenk	Hay flower	1400A120116	1
Klenk	Hay flower	A130304	1
Klenk	Heather flower, rubbed	13014202	1
Klenk	Henna powder	2074A111004	1
Klenk	Henna powder	20747A111004	1
Klenk	Hop strobile	1451B-120817	1
Klenk	Hop strobile	1451b130409	1
Klenk	Hop strobile	1451B120817	7
Klenk	Hop strobile	1451B111004	1
Klenk	Hop strobile	1450A110712	1
Klenk	Hop strobile	1451B120202	6
Klenk	Hop strobile	1451B140408	1
Klenk	Hop strobile	B130409	1
Klenk	Hop strobile	1451B130409	1
Klenk	Hop strobile	1451b130107	1
Klenk	Hop strobile	20071301	1
Klenk	Hop strobile	1451B130409	8
Klenk	Horse chestnut flower...	1440A110311	3
Klenk	Horse chestnut leaf, ...	2551B081223	1
Klenk	Horse chestnut leaf, ...	N130124	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Horse chestnut leaf, ...	255100903	1
Klenk	Horse chestnut leaf, ...	2551N130124	3
Klenk	Iceland moss, cut	5191A140828	1
Klenk	Iceland moss, cut	5191N110609	1
Klenk	Ispaghula seed	7560A140203	1
Klenk	Ispaghula seed	7560A120828	7
Klenk	Ispaghula seed	7560A130418	2
Klenk	Ispaghula seed	7560A130710	1
Klenk	Ispaghula seed	7560A150420	3
Klenk	Ispaghula seed	20071302	1
Klenk	Ispaghula seed	7560A160229	1
Klenk	Ivy leaf, cut	B130118	1
Klenk	Jasmine flower	1460A110505	1
Klenk	Java tea, cut	2221B110718	6
Klenk	Java tea, cut	17121304	1
Klenk	Java tea, cut	2221B12112	1
Klenk	Java tea, cut	B121112	3
Klenk	Java tea, cut	B110718	1
Klenk	Juniper	2890A130513	2
Klenk	Juniper	2890A101201	5
Klenk	Juniper	2890a120709	1
Klenk	Juniper	A120709	3
Klenk	Juniper	2890A120709	13
Klenk	Juniper, pounded	2896B121016	7
Klenk	Kaolin, heavy	45131I01190	1
Klenk	Kaolin, heavy	424A130507	2
Klenk	Kaolin, heavy	424A100819	2
Klenk	Kelp, cut	1051A110920	1
Klenk	Knotgrass, cut	4301A121120	3
Klenk	Knotgrass, cut	4301A120912	1
Klenk	Ladys bedstraw herb,...	3921C130226	1
Klenk	Ladys bedstraw herb,...	3921B12030	1
Klenk	Lapacho wood, cut	13298302	1
Klenk	Large waybread herb, ...	4281N110404	1
Klenk	Large waybread herb, ...	4281N111121	5
Klenk	Large waybread herb, ...	N111121	1
Klenk	Laurel leaf	B110809	1
Klenk	Laurel leaf, cut	B110809	2
Klenk	Lavender flower	1489A100125	1
Klenk	Lavender flower	1488B130604	5
Klenk	Lavender flower	1489A130214	16
Klenk	Lavender flower	313U-01462	1
Klenk	Lavender flower	1488B130111	2
Klenk	Lavender flower	1488A120803	3
Klenk	Lavender flower	1488A111201	2
Klenk	Lavender flower	1489A120103	19
Klenk	Lavender flower	1489A1203	1
Klenk	Lemon grass, cut	7721A130708	1
Klenk	Lemon grass, cut	A120906	1
Klenk	Lemon grass, cut	A130907	1
Klenk	Lemon grass, cut	7721A120906	1
Klenk	Lemon peel, cut	621A130221	1
Klenk	Lemon peel, cut	621A091016	1
Klenk	Lemon peel, cut	B130521	1
Klenk	Lemon peel, cut	A130221	3
Klenk	Lemon verbena leaf, cut	2411B110408	2
Klenk	Lime flower, cut	B130905	4
Klenk	Lime flower, cut	1611a121120	2
Klenk	Lime flower, cut	1611B121105	4
Klenk	Lime flower, cut	1611A121120	7

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Lime flower, cut	Pr.Nr.4411A-01474	1
Klenk	Lime flower, cut	1611B120827	11
Klenk	Lime flower, cut	1611B110421	2
Klenk	Lime flower, cut	1611B111011	5
Klenk	Lime flower, cut	B130307	2
Klenk	Lime flower, cut	1611B130307	4
Klenk	Lime flower, cut	1611B131220	1
Klenk	Lime flower, cut	1611B130905	5
Klenk	Lime flower, cut	1611B1010004	1
Klenk	Lime leaf, cut	2471B120502	1
Klenk	Linseed	7307A121212	1
Klenk	Linseed	7307a111107	1
Klenk	Linseed	903312240	1
Klenk	Linseed	1910212240	1
Klenk	Linseed	7307A130124	1
Klenk	Linseed	7300A111107	4
Klenk	Linseed	7307A130110	15
Klenk	Linseed	7308A111107	3
Klenk	Linseed	7300A130110	7
Klenk	Linseed	1913Q-02240	1
Klenk	Linseed	7307A160414	1
Klenk	Linseed	7308A150414	2
Klenk	Linseed	7307A150710	3
Klenk	Linseed	7300A141201	4
Klenk	Linseed	7300A140526	1
Klenk	Linseed	7500A140556	1
Klenk	Linseed	3707A140526	1
Klenk	Linseed	7307A140526	3
Klenk	Linseed	7307a131104	1
Klenk	Linseed	2002285	1
Klenk	Linseed	7307A131104	8
Klenk	Linseed	7307A131021	1
Klenk	Linseed	7307A-130722	1
Klenk	Linseed	7308A150710	2
Klenk	Linseed	7300A130722	8
Klenk	Linseed	7307A-130110	1
Klenk	Linseed	A130110	1
Klenk	Linseed	7308A130110	4
Klenk	Linseed	7300A170306	2
Klenk	Linseed	7307A130722	5
Klenk	Liquorice root, peeled	A130426	1
Klenk	Liquorice root, peeled	6381A140109	2
Klenk	Liquorice root, peeled	6381A130820	4
Klenk	Liquorice root, peeled	6381A130426	21
Klenk	Liquorice root, peeled	638A130820	1
Klenk	Liquorice root, peeled	A130820	3
Klenk	Liquorice root, peeled	6381A120806	3
Klenk	Liquorice root, peeled	6381a121126	1
Klenk	Liquorice root, peeled	6381A110531	2
Klenk	Liquorice root, peeled	6381A130425	1
Klenk	Liquorice root, peeled	6381a111004	1
Klenk	Liquorice root, peeled	6381A121126	27
Klenk	Liquorice root, natur...	6391A080905	4
Klenk	Liquorice root, natur...	6391A110421	4
Klenk	Liquorice root, natur...	6390A101020	1
Klenk	Lovage root, cut	6371B121029	3
Klenk	Lovage root, cut	B121029	1
Klenk	Lovage root, cut	B110228	1
Klenk	Lovage root, cut	6371B-121029	1
Klenk	Lovage root, cut	6371A080508	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Lungwort herb, cut	4331A120813	1
Klenk	Lungwort herb, cut	4331B111205	1
Klenk	Lycopodium herb, cut	4111B110428	1
Klenk	Lycopodium herb, cut	4111B120220	1
Klenk	Maize style, cut	A110411	1
Klenk	Maize style, cut	5581A111206	1
Klenk	Maize style, cut	5581A081204	1
Klenk	Male speedwell wort, ...	6661C130226	1
Klenk	Male speedwell wort, ...	4661c120321	1
Klenk	Male speedwell wort, ...	4661C130226	27
Klenk	Male speedwell wort, ...	4661C131104	1
Klenk	Male speedwell wort, ...	4661C2004	1
Klenk	Male speedwell wort, ...	4661C130226-32743106	1
Klenk	Male speedwell wort, ...	4661C13026	1
Klenk	Male speedwell wort, ...	4661B101011	1
Klenk	Male speedwell wort, ...	4661B110121	2
Klenk	Male speedwell wort, ...	4661B110916	1
Klenk	Male speedwell wort, ...	4661C111011	13
Klenk	Male speedwell wort, ...	C130226	5
Klenk	Male speedwell wort, ...	4661c120802	1
Klenk	Male speedwell wort, ...	4661C-111011	1
Klenk	Male speedwell wort, ...	C120321	1
Klenk	Male speedwell wort, ...	4661C120802	13
Klenk	Male speedwell wort, ...	4661A090318	1
Klenk	Male speedwell wort, ...	4661c130226-34343111	1
Klenk	Mallow flower, cut	1910211464	1
Klenk	Mallow flower, cut	1501B100408	1
Klenk	Mallow flower, cut	1501B120427	10
Klenk	Mallow flower, cut	10357908	2
Klenk	Mallow flower, cut	1501B130829	3
Klenk	Mallow flower, cut	B130829	3
Klenk	Mallow flower, cut	B120427	1
Klenk	Mallow flower, cut	150B111215	1
Klenk	Mallow flower, cut	1501b111215	1
Klenk	Mallow flower, cut	1501b110207	1
Klenk	Mallow flower, cut	1501B111215	1
Klenk	Mallow flower, cut	1501B-120427	1
Klenk	Mallow leaf, cut	B130225	1
Klenk	Mallow leaf, cut	2131B10529	1
Klenk	Mallow leaf, cut	2131B130225	4
Klenk	Manna, cut	5221A120509	10
Klenk	Manna, cut	5221A110817	1
Klenk	Manna, cut	5221A140210	4
Klenk	Manna, cut	5221A150107	2
Klenk	Mannitol	1211034-02	18
Klenk	Mannitol	5224A110314	1
Klenk	Marian thistle fruit	7180A110209	2
Klenk	Marian thistle fruit	7180A120725	1
Klenk	Marian thistle fruit	A120725	1
Klenk	Marian thistle fruit	7180A132019	1
Klenk	Marian thistle fruit	7180A130219	1
Klenk	Marian thistle fruit	7180A090929	1
Klenk	Marian thistle herb, ...	3701A111207	11
Klenk	Marian thistle herb, ...	3701A100916	2
Klenk	Marian thistle herb, ...	11041310	1
Klenk	Marian thistle herb, ...	3701B111117	1
Klenk	Marian thistle herb, ...	913E-01627	1
Klenk	Marian thistle herb, ...	2011402	1
Klenk	Marian thistle herb, ...	3701a111207	1
Klenk	Marian thistle herb, ...	3701B130318	5

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Marian thistle herb, ...	B110728	1
Klenk	Marian thistle herb, ...	3701B120717	3
Klenk	Marjoram herb, rubbed	4881A101025	15
Klenk	Marjoram herb, rubbed	4881a101025	1
Klenk	Marjoram herb, rubbed	4881A1311105	1
Klenk	Marjoram herb, rubbed	12345678	1
Klenk	Marjoram herb, rubbed	A101025	2
Klenk	Marrubium vulgare her...	13289102	1
Klenk	Marrubium vulgare her...	4141A130110	1
Klenk	Marshmallow root, pee...	6071A120314	10
Klenk	Marshmallow root, pee...	A130408	3
Klenk	Marshmallow root, pee...	6071A-130408	1
Klenk	Marshmallow root, pee...	6074B100416	1
Klenk	Marshmallow root, pee...	6071A-120314	1
Klenk	Marshmallow root, pee...	6071B091207	1
Klenk	Marshmallow root, pee...	130116	1
Klenk	Marshmallow root, pee...	6071A130612	4
Klenk	Marshmallow root, pee...	6071a130408	1
Klenk	Maté leaf, cut	2141a111213	1
Klenk	Maté leaf, cut	2131A111213	1
Klenk	Maté leaf, cut	2141A111213	16
Klenk	Maté leaf, cut	2141A-111213	1
Klenk	Matricaria flower	4100113002	2
Klenk	Matricaria flower	1290A130516	9
Klenk	Matricaria flower	2310211452	1
Klenk	Matricaria flower	1290A120403	3
Klenk	Matricaria flower	27111204	1
Klenk	Matricaria flower	1290A110607	1
Klenk	Matricaria flower	1290A120620	13
Klenk	Matricaria flower	1290A120404	6
Klenk	Matricaria flower	1290A120404	1
Klenk	Matricaria flower	1290A20404	1
Klenk	Matricaria flower	1290A110421	1
Klenk	Matricaria flower	1290A140508	1
Klenk	Matricaria flower	1290A12061	1
Klenk	Matricaria flower	1290A150618	1
Klenk	Matricaria flower	1290A120601	14
Klenk	Matricaria flower	1290A-120620	1
Klenk	Matricaria flower	A120529	1
Klenk	Matricaria flower	1290A130502	1
Klenk	Matricaria flower	1290A130502	1
Klenk	Matricaria flower	1290A130516	3
Klenk	Matricaria flower	16121301	1
Klenk	Matricaria flower	A120529	1
Klenk	Matricaria flower	A130516	1
Klenk	Matricaria flower	1290a120620	1
Klenk	Matricaria flower	14111314	1
Klenk	Matricaria flower	A130502	1
Klenk	Matricaria flower	11101302	1
Klenk	Matricaria flower	B130415	2
Klenk	Matricaria flower	1290A120529	1
Klenk	Matricaria flower	1290B130415	3
Klenk	Matricaria flower	1290A120601	2
Klenk	Matricaria flower	1290B130415	38
Klenk	Matricaria flower	1290A120529	5
Klenk	Matricaria flower	1290A110825	1
Klenk	Matricaria flower	1290B140123	2
Klenk	Matricaria flower	1290B140624	1
Klenk	Meadowsweet flower, r...	A120312	1
Klenk	Melilot herb, cut	4151A110324	8

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Melilot herb, cut	4151A120912	1
Klenk	Melilot herb, cut	A110324	6
Klenk	Melissa leaf, cut	2161A111207	15
Klenk	Melissa leaf, cut	2161A101216	1
Klenk	Melissa leaf, cut	2161A151209	1
Klenk	Melissa leaf, cut	2161A121206	2
Klenk	Melissa leaf, cut	2161A121009/PN:1198	1
Klenk	Melissa leaf, cut	2161A121029	4
Klenk	Melissa leaf, cut	2161A121009	14
Klenk	Melissa leaf, cut	2161A120718	11
Klenk	Melissa leaf, cut	27091201	1
Klenk	Melissa leaf, cut	2161A110818	1
Klenk	Melissa leaf, cut	2807211503	1
Klenk	Melissa leaf, cut	111207	1
Klenk	Mistletoe, cut	2901311695	1
Klenk	Mistletoe, cut	4691A121122-35454110	1
Klenk	Mistletoe, cut	4691a130307	1
Klenk	Mistletoe, cut	4691a121122	1
Klenk	Mistletoe, cut	4961A120419	1
Klenk	Mistletoe, cut	110816	1
Klenk	Mistletoe, cut	4691A121122	10
Klenk	Mistletoe, cut	4691A120507	2
Klenk	Mistletoe, cut	4691A130307	7
Klenk	Mistletoe, cut	1512211695	1
Klenk	Mistletoe, cut	4691A120419	6
Klenk	Mistletoe, cut	4691A-120507	1
Klenk	Mistletoe, cut	4691A110816	2
Klenk	Motherwort herb, cut	4071A-100707	1
Klenk	Mullein flower, cut	1651B150722	1
Klenk	Mullein flower, cut	1651B130820	5
Klenk	Mullein flower, cut	1650A111018	3
Klenk	Mullein flower, cut	1651B120315	3
Klenk	Mullein flower, cut	1651B140324	2
Klenk	Mullein flower, cut	1651B121109	22
Klenk	Mullein flower, cut	A120612	1
Klenk	Mullein flower, cut	1650A120612	1
Klenk	Mullein flower, cut	1650A100608	1
Klenk	Mullein flower, cut	1651b130820	2
Klenk	Mullein flower, cut	1651B110203	2
Klenk	Myrrh	3351A160330	2
Klenk	Myrrh	3351A121004	1
Klenk	Myrrh	3351A130902	1
Klenk	Myrrh	3351A120514	1
Klenk	Nettle leaf, cut	A110913	1
Klenk	Nettle leaf, cut	2381a130122	1
Klenk	Nettle leaf, cut	1107211689	1
Klenk	Nettle leaf, cut	2381A130703	13
Klenk	Nettle leaf, cut	A130703	1
Klenk	Nettle leaf, cut	4513E3054	1
Klenk	Nettle leaf, cut	A130122	1
Klenk	Nettle leaf, cut	2381A141124	1
Klenk	Nettle leaf, cut	2381A151012	1
Klenk	Nettle leaf, cut	2381A140221	1
Klenk	Nettle leaf, cut	2381A110809	1
Klenk	Nettle leaf, cut	2381A-121010	1
Klenk	Nettle leaf, cut	2381A-120619	1
Klenk	Nettle leaf, cut	2381A150109	1
Klenk	Nettle leaf, cut	2381a120619	4
Klenk	Nettle leaf, cut	1913A-03054	1
Klenk	Nettle leaf, cut	2381A121010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Nettle leaf, cut	114Q-03054	1
Klenk	Nettle root, cut	6631N110513	1
Klenk	Nettle root, cut	9931B120927	1
Klenk	Nettle seed	7470A170130	1
Klenk	Nettle seed	7470A160928	2
Klenk	Nettle seed	7470A151009	1
Klenk	Nettle seed	2021505	1
Klenk	Nettle seed	7470A110418	1
Klenk	Nutmeg, powdered (spi...	A130304	2
Klenk	Nutmeg, powdered (spi...	7325A130304	1
Klenk	Nutmeg, powdered (spi...	7325A130515	1
Klenk	Nutmeg, powdered (spi...	7325A120903	2
Klenk	Oak bark, cut	761A110816	1
Klenk	Oak bark, cut	C-280513-2	1
Klenk	Oak bark, cut	761A110816	10
Klenk	Oak bark, cut	761A110816	1
Klenk	Oak bark, cut	761A090420	1
Klenk	Olive leaf	2211B120417	2
Klenk	Olive leaf	2211B131025	1
Klenk	Olive leaf	2211B150818	1
Klenk	Olive leaf	2211B130425	3
Klenk	Orange peel, cut	531A091125	3
Klenk	Orange peel, cut	531a130204	1
Klenk	Orange peel, cut	531A130204	5
Klenk	Oregano herb	4221A120315	2
Klenk	Orris root	6831B111013	1
Klenk	Orris root	B111013	1
Klenk	Orris root	B130102	1
Klenk	Oxybuprocaine hydroch...	3613M-07321	1
Klenk	Parsley root, cut	6441A121024	4
Klenk	Parsley root, cut	6441A120283	1
Klenk	Parsley root, cut	6441A150811	1
Klenk	Parsley root, cut	6441A120606	5
Klenk	Parsley seed	7370A091120	2
Klenk	Parsley seed	7370A09120	1
Klenk	Passionflower herb, cut	4231A120329	2
Klenk	Passionflower herb, cut	A120329	2
Klenk	Passionflower herb, cut	4231A-110913	1
Klenk	Peony root, cut	B110916	2
Klenk	Peony root, cut	6421B120327	1
Klenk	Peppermint leaf, cut	2181A-120913	1
Klenk	Peppermint leaf, cut	2161A110207	1
Klenk	Peppermint leaf, cut	1913A-01505	1
Klenk	Peppermint leaf, cut	2161A121009	6
Klenk	Peppermint leaf, cut	A121009	1
Klenk	Peppermint leaf, cut	2161A120718	3
Klenk	Peppermint leaf, cut	2183A120913	2
Klenk	Peppermint leaf, cut	2161	1
Klenk	Peppermint leaf, cut	2181A130320	6
Klenk	Peppermint leaf, cut	2181A100323	1
Klenk	Peppermint leaf, cut	216A121029	1
Klenk	Peppermint leaf, cut	2181a120913	3
Klenk	Peppermint leaf, cut	2161A121029	16
Klenk	Peppermint leaf, cut	2181A120913	34
Klenk	Peppermint leaf, cut	4513E1503	1
Klenk	Peppermint leaf, cut	4113Q1505	1
Klenk	Peppermint leaf, cut	A130806	1
Klenk	Peppermint leaf, cut	2161A07020203	1
Klenk	Peppermint leaf, cut	2181B130717	10
Klenk	Peppermint leaf, cut	1001311505	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Peppermint leaf, cut	2181A101013	6
Klenk	Peppermint leaf, cut	2181A-101013	1
Klenk	Peppermint leaf, cut	1209211505	1
Klenk	Peppermint leaf, cut	CH.2161B080818	1
Klenk	Peppermint leaf, cut	2161A130806	7
Klenk	Peppermint leaf, cut	2181A13020	1
Klenk	Peppermint leaf, cut	2161A-121206	1
Klenk	Peppermint leaf, cut	2161B140321	1
Klenk	Peppermint leaf, cut	2181A151119	1
Klenk	Peppermint leaf, cut	2161A111207	1
Klenk	Peppermint leaf, cut	2161A121206	14
Klenk	Potato starch	94A110512	1
Klenk	Potato starch	94A130304	4
Klenk	Potato starch	94a170123	1
Klenk	Potato starch	94A170123	2
Klenk	Primrose flower	1541B110314	7
Klenk	Primrose flower	1540A100907	1
Klenk	Primrose flower	B130808	1
Klenk	Primrose flower	1541B100414	1
Klenk	Primrose flower	1540A-100907	1
Klenk	Primrose flower	A111018	1
Klenk	Primrose flower	B130808	1
Klenk	Primrose flower	B110314	1
Klenk	Primrose root, cut	B130218	3
Klenk	Primrose root, cut	6461A100511	1
Klenk	Progesterone, microni...	C01-120322	2
Klenk	Psyllium seed	7380a120109	3
Klenk	Psyllium seed	7380A150520	4
Klenk	Psyllium seed	7380A130128	6
Klenk	Psyllium seed	A130128	1
Klenk	Psyllium seed	A120109	1
Klenk	Psyllium seed	7380A120911	16
Klenk	Psyllium seed	7380A110317	1
Klenk	Psyllium seed	7560a120330	1
Klenk	Psyllium seed	7380A111018	1
Klenk	Psyllium seed	7380A120109	24
Klenk	Psyllium seed	7380a111018	1
Klenk	Psyllium seed	7560A140203	1
Klenk	Psyllium seed	A120911	1
Klenk	Psyllium seed	7560A130418	1
Klenk	Ramsons herb, cut	3511B111005	2
Klenk	Raspberry leaf, cut	A130430	1
Klenk	Raspberry leaf, cut	2291A110125	1
Klenk	Raspberry leaf, cut	2291A121120	4
Klenk	Raspberry leaf, cut	2291A120813	2
Klenk	Raspberry leaf, cut	2291A130320	1
Klenk	Raspberry leaf, cut	12273107	1
Klenk	Raspberry leaf, cut	2291A-120806	1
Klenk	Raspberry leaf, cut	2291A120507	5
Klenk	Raspberry leaf, cut	2291a111124	1
Klenk	Raspberry leaf, cut	2291A111124	1
Klenk	Red clover flower	1640A110727	1
Klenk	Red clover flower	7470A130311	1
Klenk	Red clover flower	1640a130122	1
Klenk	Red clover flower	1640A130122	15
Klenk	Red clover flower	811Q-07163	1
Klenk	Red-rose petal	1561A121016	1
Klenk	Red-rose petal, cut	1561A110404	1
Klenk	Red-rose petal, cut	1561A121016-34023101	1
Klenk	Red-rose petal, cut	1561B101210	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Redwood	5961B111219	1
Klenk	Restharrow root, cut	6411111110	1
Klenk	Restharrow root, cut	6411B101216	2
Klenk	Restharrow root, cut	6411B130404	2
Klenk	Restharrow root, cut	4296345	1
Klenk	Restharrow root, cut	6411B111110	1
Klenk	Restharrow root, cut	B111110	1
Klenk	Rhubarb, cut	6891B111118	1
Klenk	Rhubarb, cut	7761B130514	2
Klenk	Ribwort plantain, cut	A110816	1
Klenk	Ribwort plantain, cut	4271A110816	2
Klenk	Ribwort plantain, cut	4271A-110816	1
Klenk	Ribwort plantain, cut	4271A111125	2
Klenk	Ribwort plantain, cut	4271A100615	2
Klenk	Ribwort plantain, cut	4281N111121	1
Klenk	Ribwort plantain, cut	4271A130716	1
Klenk	Ribwort plantain, cut	4271A130115	6
Klenk	Ribwort plantain, cut	427A111125	1
Klenk	Rice starch	84A130304	1
Klenk	Rice starch	64a130304	1
Klenk	Rice starch	94A100122	1
Klenk	Rice starch	104A120621	5
Klenk	Rice starch	64A080218	1
Klenk	Rice starch	104A141114	2
Klenk	Rice starch	104A130304	8
Klenk	Rice starch	64A130304	11
Klenk	Rice starch	84A141114	4
Klenk	Robert herb, cut	3951B12118	1
Klenk	Robert herb, cut	3951A110621	7
Klenk	Robert herb, cut	3951B121218	2
Klenk	Robert herb, cut	A110621	6
Klenk	Rose hips peel, cut	2851A-110202	1
Klenk	Rose hips peel, cut	2851A101119	3
Klenk	Rose hips peel, cut	2851A120102	5
Klenk	Rose hips peel, cut	2851A120904	2
Klenk	Rose hips peel, cut	3613M-01543	1
Klenk	Rose hips peel, cut	2851A150310	4
Klenk	Rose hips peel, cut	A120904	3
Klenk	Rose hips peel, cut	2851A110202	1
Klenk	Rose hips peel, cut	2841B111019	2
Klenk	Rose hips with peel, ...	A120719	1
Klenk	Rose hips with peel, ...	3613M-01542	1
Klenk	Rose hips with peel, ...	A120904	1
Klenk	Rose hips with peel, ...	2841B081210	1
Klenk	Rose hips with peel, ...	2841B111019	3
Klenk	Rose hips with peel, ...	B130523	1
Klenk	Rose hips with peel, ...	2851A120904	1
Klenk	Rose hips with peel, ...	B111019	1
Klenk	Roselle	1431B130205	1
Klenk	Roselle	1431B131015	2
Klenk	Roselle	B130205	1
Klenk	Roselle	1431B-110519	1
Klenk	Roselle	1431B111122	1
Klenk	Rosemary leaf, cut	2261A160905	2
Klenk	Rosemary leaf, cut	2261A081117	1
Klenk	Rosemary leaf, cut	112E-01510	1
Klenk	Rosemary leaf, cut	301401	1
Klenk	Rosemary leaf, cut	2261A130201	29
Klenk	Rosemary leaf, cut	2261A120224	18
Klenk	Rosemary leaf, cut	2261a160329	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Rosemary leaf, cut	2261A150622	6
Klenk	Rosemary leaf, cut	A130201	2
Klenk	Rosemary leaf, cut	2261A131118	12
Klenk	Rosemary leaf, cut	2261A160329	2
Klenk	Rosemary leaf, cut	2261A140910	4
Klenk	Rosemary leaf, cut		1
Klenk	Rosemary leaf, cut	2261a120224	2
Klenk	Rosemary leaf, cut	61A131118	1
Klenk	Rue herb, cut	4361A100707	3
Klenk	Rue herb, cut	4361A07020601	1
Klenk	Rue herb, cut	4361A120207	1
Klenk	Rue herb, cut	A120207	1
Klenk	Saffron (spice)	384A-120410	1
Klenk	Saffron (spice)	384a120410	1
Klenk	Saffron (spice)	384A110608	1
Klenk	Saffron (spice)	384A130320	2
Klenk	Saffron filaments	380A130613	1
Klenk	Sage leaf, cut	2301A110420	1
Klenk	Sage leaf, cut	2301a120403	1
Klenk	Sage leaf, cut	C120126	1
Klenk	Sage leaf, cut	2301C120126	16
Klenk	Sage leaf, cut	2301C12026	1
Klenk	Sage leaf, cut	2301c120126	1
Klenk	Sage leaf, cut	2301A110420	2
Klenk	Sage leaf, cut	2301C12126	1
Klenk	Salicylic acid	12134722	5
Klenk	Sarsaparilla root, cut	B111108	1
Klenk	Sarsaparilla root, cut	6561B111108	1
Klenk	Senna leaf, cut	B130214	1
Klenk	Senna leaf, cut	2331B130214	5
Klenk	Senna leaf, cut	2331C120614	8
Klenk	Senna leaf, cut	2331A100525	2
Klenk	Senna leaf, cut	2331A110207	1
Klenk	Senna leaf, cut	2331B120403	1
Klenk	Senna leaf, cut	2331C130326	16
Klenk	Senna leaf, cut	2331c130326	1
Klenk	Senna leaf, cut	2331B110907	2
Klenk	Senna leaf, cut	2331B150903	1
Klenk	Senna leaf, powdered	2331B150903	1
Klenk	Senna pods, tinnevelly	1041B120327	2
Klenk	Shepherd's purse, cut	3671A-100414	2
Klenk	Shepherd's purse, cut	3671B111201	1
Klenk	Shepherd's purse, cut	B111201	1
Klenk	Silverweed, cut	3551A130903	1
Klenk	Silverweed, cut	H-280513-3	1
Klenk	Silverweed, cut	11328103	4
Klenk	Silverweed, cut	A110913	2
Klenk	Solanum dulcamara ste...	5601A080513	3
Klenk	Spray-dried arabic gum	3314A140121	1
Klenk	St. John's wort, cut	4021a120604	1
Klenk	St. John's wort, cut	4021A00819	1
Klenk	St. John's wort, cut	4021B110321	2
Klenk	St. John's wort, cut	A130612	3
Klenk	St. John's wort, cut	4021A100919	1
Klenk	St. John's wort, cut	4021A120604	15
Klenk	St. John's wort, cut	15068103	2
Klenk	St. John's wort, cut	4021A100819	1
Klenk	St. John's wort, cut	A121120	1
Klenk	St. John's wort, cut	4021A120905	8
Klenk	St. John's wort, cut	4021A130612	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	St. John's wort, cut	4021A-100819	1
Klenk	St. John's wort, cut	4021A100819	21
Klenk	St. John's wort, cut	4021A100819914	2
Klenk	St. John's wort, cut	4021a100819	1
Klenk	St. John's wort, cut	4012A100819	1
Klenk	Star anise	2690A111013	3
Klenk	Star anise	2690A120809	1
Klenk	Star anise	B111109	1
Klenk	Star anise	2690A101026	1
Klenk	Star anise	2696B111109	1
Klenk	Star anise	2690A130903	2
Klenk	Stemless carline root	69191B120726	1
Klenk	Stemless carline root	B110831	1
Klenk	Stemless carline root	6194b110831	1
Klenk	Stemless carline root	12097301	1
Klenk	Stemless carline root	6191B120726	3
Klenk	Stemless carline root	6194B110831	1
Klenk	Strawberry leaf, cut	A130430	1
Klenk	Strawberry leaf, cut	2031A130430	2
Klenk	Strawberry leaf, cut	B120925	1
Klenk	Strawberry leaf, cut	2031B110224	2
Klenk	Sunflower, cut	1711I-07036	1
Klenk	Sunflower, cut	1701B091013	1
Klenk	Sweet fennel	2870A131002	1
Klenk	Sweet fennel	2870A140528	2
Klenk	Sweet fennel	2870130304A	1
Klenk	Sweet woodruff herb, ...	3601B140321	1
Klenk	Sweet woodruff herb, ...	3601B090114	1
Klenk	Sweet woodruff herb, ...	B090114	3
Klenk	Tall golden rod, cut	4521A110209	4
Klenk	Tall golden rod, cut	4521A-110816	1
Klenk	Tall golden rod, cut	A110209	1
Klenk	Tall golden rod, cut	4521A100414	1
Klenk	Tall hawkweed, cut	B111005	1
Klenk	Tall hawkweed, cut	3611B11105	1
Klenk	Tartaric acid	13405	2
Klenk	Theriac, powdered	12189403	1
Klenk	Thyme, rubbed	4622B110519	1
Klenk	Thyme, rubbed	4622a120111	3
Klenk	Thyme, rubbed	1910211688	1
Klenk	Thyme, rubbed	4622A141203	1
Klenk	Thyme, rubbed	4622A110728	25
Klenk	Thyme, rubbed	4622A120111	35
Klenk	Thyme, rubbed	4622a110728	1
Klenk	Thyme, rubbed	4622A-120111	1
Klenk	Thyme, rubbed	A130910	3
Klenk	Thyme, rubbed	4622A130910	2
Klenk	Thyme, rubbed	4622a130109	2
Klenk	Thyme, rubbed	2103311688	1
Klenk	Tormentil	6911B111013	10
Klenk	Tormentil	6911B140227	1
Klenk	Tormentil	6911B141114	2
Klenk	Tormentil	6914b120112	1
Klenk	Tormentil	6911B150825	4
Klenk	Tormentil	6911B131029	4
Klenk	Turmeric root, cut	6741B121030	1
Klenk	Turmeric, javanese, cut	B130117	1
Klenk	Turmeric, javanese, cut	6731B150630	1
Klenk	Valerian root, cut	6645B130226	1
Klenk	Valerian root, cut	13196808	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Valerian root, cut	6641B-110317	1
Klenk	Valerian root, cut	6641B130226	9
Klenk	Valerian root, cut	6641B110317	2
Klenk	Valerian root, cut	6641B100823	1
Klenk	Valerian root, cut	6641B120417	1
Klenk	Verbena herb, cut	4651A070828	3
Klenk	Verbena herb, cut	4651B130611	2
Klenk	Verbena herb, cut	A130923	1
Klenk	Verbena herb, cut	4654B120611	1
Klenk	Verbena herb, cut	4651B120611	17
Klenk	Walnut bark, cut	711B080806	1
Klenk	Walnut bark, cut	711B120613	3
Klenk	Walnut bark, cut	711b120613	1
Klenk	Walnut bark, cut	711B120906	1
Klenk	Walnut leaf, cut	A130410	1
Klenk	Walnut leaf, cut	2111A-120123	1
Klenk	Walnut leaf, cut	F-280513-1	1
Klenk	Walnut leaf, cut	2111A110816	1
Klenk	Walnut leaf, cut	2111A131029	1
Klenk	Walnut leaf, cut	2111A130410	2
Klenk	Wax, bleached	280A160311	12
Klenk	Wax, bleached	280A150528	7
Klenk	Wax, bleached	290A150716	1
Klenk	Wax, bleached	280A171019	11
Klenk	Wax, yellow	13335513	4
Klenk	Wax, yellow	290A151214	3
Klenk	Wax, yellow	280A110110	1
Klenk	Wax, yellow	290A120627	5
Klenk	Wax, yellow	290A110803	1
Klenk	Wax, yellow	282A120627	1
Klenk	Wax, yellow	280-A120627	1
Klenk	Wax, yellow	290A150716	3
Klenk	Wax, yellow	282A140306	1
Klenk	Wax, yellow	280A150528	5
Klenk	Wax, yellow	280A130115	28
Klenk	Wax, yellow	290A081210	1
Klenk	Wax, yellow	280A140306	8
Klenk	Wax, yellow	290A140903	5
Klenk	Wax, yellow	280A13011A	1
Klenk	Wax, yellow	160311	1
Klenk	Wax, yellow	280A171019	1
Klenk	Wax, yellow	290A170915	5
Klenk	Wax, yellow	290A170314	3
Klenk	Wax, yellow	190A170314	1
Klenk	Wheat starch	104a141114	1
Klenk	Wheat starch	104A141114	10
Klenk	Wheat starch	104A101116	1
Klenk	Wheat starch	104A160512	9
Klenk	Wheat starch	104A171122	2
Klenk	Wheat starch	10A160512	1
Klenk	White deadnettle flower	1470A080103	1
Klenk	White deadnettle flower	1470A120111	4
Klenk	White mustard seed, p. . .	7254A141202	1
Klenk	Wild pansy, cut	4681B130715	5
Klenk	Wild pansy, cut	4681A-120806	1
Klenk	Wild pansy, cut	4681A120806	1
Klenk	Wild pansy, cut	4681A110303	2
Klenk	Wild pansy, cut	4681A120806	15
Klenk	Wild pansy, cut	B130715	1
Klenk	Willow bark, cut	811B100910	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Klenk	Willow bark, cut	B100910	6
Klenk	Willowherb, small flo...	4721C2120321	1
Klenk	Willowherb, small flo...	4721A120418	11
Klenk	Willowherb, small flo...	4721C110701	4
Klenk	Willowherb, small flo...	2512U-01641	1
Klenk	Willowherb, small flo...	A120418	1
Klenk	Willowherb, small flo...	4721C120321	10
Klenk	Willowherb, small flo...	4721C130618	4
Klenk	Willowherb, small flo...	4721A101010	2
Klenk	Willowherb, small flo...	4721A120903	7
Klenk	Willowherb, small flo...	4721B130528	1
Klenk	Willowherb, small flo...	C120321	1
Klenk	Wormwood herb, cut	3461A1100324	1
Klenk	Wormwood herb, cut	A110324	1
Klenk	Wormwood herb, cut	3461A121130	2
Klenk	Wormwood herb, cut	3461A110324	2
Klenk	Yarrow flower, rubbed	4181A121130	1
Klenk	Yarrow flower, rubbed	A110816	3
Klenk	Yarrow herb, cut	4181A121130	12
Klenk	Yarrow herb, cut	4181A101102	1
Klenk	Yarrow herb, cut	4181a130513	1
Klenk	Yarrow herb, cut	4181A130710	1
Klenk	Yarrow herb, cut	A130709	1
Klenk	Yarrow herb, cut	4181B130709	2
Klenk	Yarrow herb, cut	4181a11011	1
Klenk	Yarrow herb, cut	15333007	1
Klenk	Yarrow herb, cut	4181130709	1
Klenk	Yarrow herb, cut	4181A110111	15
Klenk	Yellow chaste weed, cut	1591A-110511	1
Klenk	Yellow chaste weed, cut	1591A110511	13
Klenk	Zedoary root, cut	6941B121025	3
Klenk	Zedoary root, cut	6941B110329	2
Klenk	Zedoary root, powdered	12007505	1
Klosterring	Betamethasone dipropi...	17B16-B04-332219	2
Klosterring	Clotrimazole	6234505	1
Klosterring	Gentamicin sulfate	15G29-B09	1
Klosterring	Metoprolol tartrate	17G03-B09-345301	1
Klosterring	Triclosan	18081405	16
Klosterring	Vanillin	17117207	7
kögl	Biotin	17315206	1
kögl	Lidocaine hydrochlori...	1803832	1
kögl	Magnesium hydrogen ci...	1900910	1
Kögl	Estradiol hemihydrate	296109164216	1
Kögl	Mannitol	1900820	1
Kögl	Prasterone	1803108	1
KÖGL	Ascorbic acid	1900900	1
KÖGL	Biotin	17315212	1
KÖGL	Macrogols 4000	1900232	1
KÖGL	Quinine sulfate dihyd...	1901102	1
Koisser	Cellulose, microcryst...	E0790	1
Koisser	Melatonin	140202	1
Koisser	Melatonin	140202	1
Koisser	Pregnenolone	ZH-130802	1
Koisser	Procaine hydrochloride	11370218	1
Koisser	Rice starch	WE07032013	1
Koisser Pharma	Pregnenolone	ZH-131001	1
KOKO Gmbh	Zinc gluconate	1810014-01	28
Kottas	Bitter fennel	P15301479	1
kräuter schulte	Bitter orange peel, cut	51265	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
kräuter schulte	Gentian root, cut	52348	1
kräuter schulte	Red-rose petal, cut	47543	1
kräuter schulte	Tormentil	51952	1
Kräuter Schulte	Black cohosh, cut	3308M-02176	1
Kräuter Schulte	Blackthorn flower, ru...	53233	2
Kräuter Schulte	Bucho leaf	54519	1
Kräuter Schulte	China Pai Mu Tan	52783	1
Kräuter Schulte	Couch grass rhizome, ...	53016	2
Kräuter Schulte	Dill fruit	6363	1
Kräuter Schulte	Dryas octopetala herb...	38970	1
Kräuter Schulte	Geum urbanum herb, cut	51737	1
Kräuter Schulte	Meadowsweet flower, r...	53473	1
Kräuter Schulte	Meadowsweet flower, r...	51369	1
Kräuter Schulte	Red clover flower	53778	1
Kräuter Schulte	Red tea	52029	1
Kräuter Schulte	Saponaria root, cut	54977	1
Kräuter Schulte	Saponaria root, cut	54500	1
Kräuter Schulte	Saponaria root, cut	54524	2
Kräuter-Schulte	Peppermint leaf, cut	56254	1
Kräuter-Schulte	Sorrel herb, cut	55272	1
Kräuter-Schulte	Sweet woodruff herb, ...	52878	1
Kremer/Phönix	Triclosan	15/157/T	1
Krieger	Chlortetracycline hyd...	12C16-No1	1
Krieger	Citric acid monohydrate	613001-1-P	1
Krieger	Emulsifying cetostear...	1711001-01	1
Krieger	Erythromycin	16G06-B03-329150	4
Krieger	Erythromycin	16J24-B10-341876	1
Krieger	Macrogols 1500	1712002	1
Krieger	Salicylic acid	16050012	13
Krieger	Salicylic acid	13032017-MW-A	1
Krieger	Tretinoin	1612010-02	14
Krieger	Triclosan	12g17-N04	1
KS Pharma	Melatonin	XPS130428	1
KS Pharma	Pregnenolone	ZH-111119	1
Kuerzel	Metronidazole	1307002-02	4
Kwizda	Alfatradiol	529068	1
L10	Citric acid	9494185	1
L10	Doxylamine hydrogen s...	16F06-B02-341879	2
L10	Doxylamine hydrogen s...	16F06-B02-351516	1
L10	Doxylamine hydrogen s...	16F06B02338805	1
L10	Hydrocortisone	16H01-B01-326071	5
L10	Metronidazole	1807039-01	12
L10	Sorbitol	16353108	2
L10	Wax, yellow	17207403	8
L11	Non-ionic emulsifying...	17043701	1
L12	Ascorbyl palmitate	18B12-B09	1
L12	Betamethasone valerate	17C29-B04-334329	2
L12	Calcium citrate tetra...	1606001-03	1
L12	Calcium gluconate	15435627	1
L12	Camphor	18108901	1
L12	Cetyl palmitate 15	18149904	10
L12	Clindamycin hydrochlo...	18009403	4
L12	Clindamycin hydrochlo...	18009402	2
L12	Cocoa butter	16316508	21
L12	Doxylamine hydrogen s...	16F06-B02-338805	5
L12	Doxylamine hydrogen s...	15K04-B01-327357	8
L12	Doxylamine hydrogen s...	15K04-B01-332765	3
L12	Doxylamine hydrogen s...	16F06-B02-348870	4
L12	Eosin Y	17K21-T07-054710	1
L12	Erythromycin	18096203	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
L12	Glycerol monostearate...	16233802	2
L12	Hydroxyethylcellulose...	18022505	2
L12	Hypromellose 4000	17115008	5
L12	Isosorbide dinitrate...	15I10-B01-331286	2
L12	Paraffin, microcrysta...	17137502	2
L12	Potassium citrate	15435709	3
L12	Propranolol hydrochlo...	18D23-B09-355583	6
L12	Resocinol	16153306	1
L12	Sorbic acid	18287602	1
L12	Sorbitol	17140904	4
L12	Thymol	15172002	1
L12	Triamcinolone, micron...	17359001	1
L12	Triclosan	1628081917181722	1
L30	Clindamycin hydrochlo...	17209301	2
L30	Dequalinium chloride	1804009-01	3
L30	Doxylamine hydrogen s...	16F06-B02-343936	5
L30	Doxylamine hydrogen s...	16F06-b02-334258	1
L30	Erythromycin	17088208	1
L30	Lactose monohydrate	1703019-02	4
L30	Progesterone, microni...	17060910	1
L30	Salicylic acid	17301507	8
L31	Heparin sodium	16I26-B02-337593	1
L31	Non-ionic emulsifying...	16223104	1
Laborchemie APOLDA GmbH	Cannabidiol	APCK312001	1
Laborchemie Apolda/Je...	Urea	12345	1
Laborchemie Apolda/Ph...	Oxytetracycline hydro...	1506037-02	32
Land/Phönix	Neomycin sulfate	13062905	4
LaVita	Lidocaine	15H27-B07-318310	6
Lihn/Reformhaus	Fructose	713E-03394	2
Lihua Pharm	Prednisolone acetate	K07A20140302M	1
Lihua/phönix	Hydrocortisone acetate	K06A20140801M	1
lohmann	Calcium citrate tetra...	3045397	1
lohmann	Magnesium hydrogen ci...	3054132	1
lohmann	Potassium citrate	1051470	1
Löns Apotheke	Urea	40003	1
Löns Apotheke	Urea	20349	1
Löwen City-Apotheke B...	Doxycycline hyclate	10022014	1
Löwen City-Apotheke B...	Doxycycline hyclate	2042015	1
Löwen City-Apotheke B...	Octenidine dihydrochl...	2042015	1
Löwen City-Apotheke B...	Propranolol hydrochlo...	4012013	1
Lysoform	Tosylchloramide sodium	K-09-2019	1
Mahendra Chem./Noweda	Lidocaine hydrochlori...	1505023	1
Martinus Apotheke	Bismuth gallath, basic	keine	1
Martinus-Apotheke	Triamcinolone acetone...	13379905	16
Mediherb	Zhu Ling		1
Mehta API	Chloramphenicol	181113082	1
Merck	Ascorbyl palmitate	16130-B03	1
Merck	Ascorbyl palmitate	17E18-B06-334914	2
Merck	Caffeine	MerckReagenziensatz	1
Merck	Doxycycline hyclate	22022019	1
Merck	Heavy alkaline bismut...	75419201	1
Merck	Histamine dihydrochlo...	WXBC3407V	1
Merck	Histamine dihydrochlo...	Wxbc3407V	1
Merck	Trichloroacetic acid	12821	1
Merck	Trometamol	8386F017	1
Midas	Clobetasol propionate	17096705	14
Midas	Linseed	7308A161208	1
Mittelwand bio 2017	Wax, yellow	5022017	1
Moehs BCN S.L. / Noweda	Lidocaine	15332305	4
Moehs Catalana, S.L./...	Lidocaine hydrochlori...		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
n	Betamethasone dipropi...	13/17-N05	11
Nachprüfung	Abrotani herb, cut	3451A070919	1
Nachprüfung	Aniseed	5010U-01530	1
Nachprüfung	Black cohosh, cut	10281113	1
Nachprüfung	Blackcurrant leaf, cut	1003	1
Nachprüfung	Boric acid		1
Nachprüfung	Butylated hydroxytolu...	123456	1
Nachprüfung	Camphor	1616E-01221	1
Nachprüfung	Cetyl palmitate 15	4917I-01242	1
Nachprüfung	Cetyl palmitate 15	3012U01241	1
Nachprüfung	Clioquinol	3021U-02562	1
Nachprüfung	Clioquinol	13L26-N15	2
Nachprüfung	Estradiol benzoate	123456	1
Nachprüfung	Fenugreek, powdered	2412E02239	1
Nachprüfung	Gentamicin sulfate	123456	1
Nachprüfung	Ginkgo leaf, cut	2461A100429	10
Nachprüfung	Glucose		1
Nachprüfung	Hay flower	1911E01459	1
Nachprüfung	Heather, cut	5110A03224	1
Nachprüfung	Hydroxyethylcellulose...	13d26-n16	1
Nachprüfung	Hydroxyethylcellulose...	123456	1
Nachprüfung	Lactose monohydrate	3712Q-02210	1
Nachprüfung	Macrogols 4000	14A03-B01-289605	1
Nachprüfung	Menthol	123456	1
Nachprüfung	Methoxsalen	4211Q-03262	1
Nachprüfung	Nystatin	123456	1
Nachprüfung	Nystatin	1118E-01941	1
Nachprüfung	Oxalic acid		1
Nachprüfung	Potassium citrate	111112	1
Nachprüfung	Prasterone	120022	1
Nachprüfung	Quinolinol sulfate po...	13071104	1
Nachprüfung	Quinolinol sulfate po...	2812U03423	1
Nachprüfung	Rice starch	3709E01116	1
Nachprüfung	Sodium citrate	1205026-02	1
Nachprüfung	Sorbitan stearate		1
Nachprüfung	Stemless carline root	2812U02133	1
Nachprüfung	Tannin	5210A-01052	1
Nachprüfung	Tretinoin		1
Nachprüfung	Wild pansy, cut	10291007	1
Nachprüfung	Xylometazoline hydroc...	4317E-03337	1
Nachprüfung 2312	Hydrocortisone	10282201	1
Nachprüfung 2574	Chlortetracycline hyd...	1108051-01	2
Nachprüfung 2609	Hypromellose 4000	11H26-N01	1
Nachprüfung 2797	Hydroxyethylcellulose...	12249202	1
Nachtestung	Carbomers 35000 / 50000	alt:115/14	1
Nachtestung	Sodium benzoate	300/18	1
Nachtestung-Verlänger...	Camphor	242/17	1
Nachtestung-Verlänger...	Citric acid	278/18	1
Nachtestung-Verlänger...	Macrogols 4000	17/17	1
Nanjing Baijingyu Pha...	Chloramphenicol	13080908	6
Newchem	Triamcinolone acetonide	1303001-02	3
NEWCHEM S.P.A.	Triamcinolone acetonide	NV0994M	1
Newchem S.p.A./Caesar...	Triamcinolone acetonide	1543015	1
Newchem Spa (site)	Triamcinolone acetonide	16K07-B06-332102	5
Norbrook/Gehe	Oxytetracycline hydro...	140701203	3
Norex Flavours Pvt. L...	Menthol	15286207	1
Norex Flavours Pvt. L...	Menthol	15286205	3
Northeast Gen	Chloramphenicol	DY122127	1
Nowed 18.12.2014	Triclosan	14121210	11
noweda	Dandelion root, cut		1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
noweda	Iceland moss, cut		1
noweda	Triamcinolone acetonide	13256604	1
noweda	Walnut leaf, cut		1
Noweda	Abrotani herb, cut	3451A120312	19
Noweda	Acetylsalicylic acid	10040403	1
Noweda	Agrimony, cut	10325503	2
Noweda	Agrimony, cut	9300106	1
Noweda	Alchemilla, cut	12153302	13
Noweda	Angelica root	272720	2
Noweda	Atropine sulfate	12250403	1
Noweda	Bentonite	11304415	1
Noweda	Betamethasone dipropi...	13236803	2
Noweda	Betamethasone dipropi...	BDP014/1115	1
Noweda	Betamethasone dipropi...	13236802	5
Noweda	Betamethasone valerate	12229832	1
Noweda	Betamethasone valerate	11B21-N02	1
Noweda	Betamethasone valerate	15075812	4
Noweda	Betamethasone valerate	11131902	1
Noweda	Betamethasone valerate	BV01090117	1
Noweda	Betamethasone valerate	14H28B01302863	2
Noweda	Betamethasone, micron...	M-X01-120301m1	1
Noweda	Betamethasone, micron...	12B20-N05	4
Noweda	Budesonide, micronized	16143110	1
Noweda	Butylated hydroxytolu...	12158201	4
Noweda	Caffeine	1404026-01	1
Noweda	Caffeine citrate	CCC-1101	1
Noweda	Calcium pantothenate	14208012	1
Noweda	Calendula flower, cut	12221709	9
Noweda	Camphor	1202021-01	3
Noweda	Caraway fruit	12226605	3
Noweda	Caraway fruit	12226503	6
Noweda	Carmellose sodium (40...	12354401	1
Noweda	Chloramphenicol	14346910	13
Noweda	Chloramphenicol	11370712	5
Noweda	Chloramphenicol	1137014	1
Noweda	Chlorhexidine diacetate	17185805	2
Noweda	Chlorhexidine diacetate	1785803	1
Noweda	Chlorocresol	12111017	4
Noweda	Chlortetracycline hyd...	12c16-n01	1
Noweda	Ciclopirox olamine	12B22-N03	2
Noweda	Ciclopirox olamine	12040911	4
Noweda	Citric acid	L1202001-02	2
Noweda	Clioquinol	14359103	11
Noweda	Clioquinol	10057206	1
Noweda	Clobetasol propionate	16106720	8
Noweda	Clotrimazole	17111608	10
Noweda	Cocoa butter	11117813	2
Noweda	Coltsfoot leaf, cut	2021B111111	3
Noweda	Damiana leaf, cut	200151	1
Noweda	Dequalinium chloride	12b14-n05	1
Noweda	Dexamethasone	12033307	6
Noweda	Diclofenac sodium	11G12-N06	1
Noweda	Diclofenac sodium	15278504	1
Noweda	Diclofenac sodium	18039202	3
Noweda	Diltiazem hydrochloride	12133010	15
Noweda	Dithranol	APDT/375/2012	2
Noweda	Elecampane rhizome, cut	6341B111118	3
Noweda	Emulsifying cetostear...	13065576	1
Noweda	Equisetum stem, cut	12231401	6
Noweda	Erythromycin	11195413	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Erythromycin	12152201	10
Noweda	Erythromycin	16079606	6
Noweda	Erythromycin	11195411	2
Noweda	Erythromycin	12152108	6
Noweda	Estradiol benzoate	13277603	3
Noweda	Estradiol hemihydrate	12268008	1
Noweda	Estradiol hemihydrate	15038004	7
Noweda	Estradiol valerate	12A20-N10	4
Noweda	Estriol	11J13-No7	1
Noweda	Estriol	12032504	3
Noweda	Fructose	1205022-01	2
Noweda	Gentamicin sulfate	1305005	1
Noweda	Glucose monohydrate	1204044-01	7
Noweda	Glucose monohydrate	1210034	1
Noweda	Glycerol monostearate...	11079301	1
Noweda	Hard fat in pastilles	11345909	1
Noweda	Hard fat W35	16027801	1
Noweda	Heavy alkaline bismut...	13/01-N02	1
Noweda	Hop strobile	274725	1
Noweda	Hydrocortisone	12104303	6
Noweda	Hydrocortisone acetate	1112031-02	8
Noweda	Hydroxyethylcellulose...	1309020-03	1
Noweda	Hypromellose 4000	1315505	1
Noweda	Iodoform	10320402	1
Noweda	Iodoform	12070407	2
Noweda	Ketamine hydrochloride	14/03-B02-311706	1
Noweda	Ketamine hydrochloride	292414P	1
Noweda	Lidocaine hydrochlori...	1812015-04	1
Noweda	Lidocaine hydrochlori...	16222207	6
Noweda	Linseed	7307A111107	7
Noweda	Liquorice root, peeled	70226/1	1
Noweda	Macrogols 4000	12361502	5
Noweda	Macrogols 4000	13136003	4
Noweda	Macrogols 4000	12361512	4
Noweda	Magnesium hydrogen ci...	10163306	1
Noweda	Manna, powdered	11170102	1
Noweda	Mannitol	P101P102P103	3
Noweda	Mannitol	11219101	2
Noweda	Mannitol	16L15-B05-329316	1
Noweda	Mannitol	16D20-B04-320487	2
Noweda	Mannitol	11219138	1
Noweda	Mannitol	15B05-B06-304302	1
Noweda	Mannitol	12105404	4
Noweda	Mannitol	15I24-B05-314323	1
Noweda	Mannitol	15I24-B05-314323	1
Noweda	Mannitol	15K17-B01	1
Noweda	Marjoram herb, rubbed	4881A100929	3
Noweda	Marrubium vulgare her...	4141A080513	1
Noweda	Menthol	15363809	9
Noweda	Menthol	17188904	10
Noweda	Methoxsalen	11D07-N16	3
Noweda	Methylrosanilinium ch...	12H28-N06	1
Noweda	Methylrosanilinium ch...	2159	2
Noweda	Metronidazole	1619-B04-329166	1
Noweda	Metronidazole	12G04-N15	18
Noweda	Metronidazole	17058005	13
Noweda	Metronidazole	160902802	1
Noweda	Metronidazole	12110404	1
Noweda	Mistletoe, cut	11294001	1
Noweda	Molecular distilled g...	11160010	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Nettle leaf, cut	12195001	7
Noweda	Nicotinamide	12133403	3
Noweda	Oak bark, cut	10348617	8
Noweda	Oat straw, cut	11280003	1
Noweda	Oxytetracycline hydro...	140701202	1
Noweda	Paraffin, hard	11363206	3
Noweda	Paraffin, microcrysta...	12266304	1
Noweda	Peppermint leaf, cut	2911U-01505	1
Noweda	Phenylephrine hydroch...	12F01-N20	2
Noweda	Potassium hydrocarbon...	1004019-01	1
Noweda	Potassium sorbate	13065222	2
Noweda	Povidone iodine	G0354BPZ0	1
Noweda	Prednisolone acetate	K07D20150105M	3
Noweda	Prednisolone acetate	1523-B04-334688	1
Noweda	Prednisolone acetate	17010410203707556	1
Noweda	Prednisolone, microni...	1108055-05	1
Noweda	Primrose root, cut	12064504	4
Noweda	Procaine hydrochloride	1309N08	1
Noweda	Propranolol hydrochlo...	155.288	1
Noweda	Quinolinol sulfate po...	140804	2
Noweda	Resocrinol	1318b09291533	1
Noweda	Resocrinol	14G10-B06-296	1
Noweda	Rice starch	14023203	4
Noweda	Rosemary leaf, cut	2261A100916	1
Noweda	Rosemary leaf, cut	15428011	1
Noweda	Salbutamol sulfate	13K12-N02	1
Noweda	Salicylic acid	12036226	8
Noweda	Sodium citrate	14A13-B02-294541	1
Noweda	Sodium citrate	11127001	1
Noweda	Sodium edetate	16337312	15
Noweda	Spironolactone	1602006-01	2
Noweda	Spironolactone	1101013-0	1
Noweda	St. John's wort, cut	12092202	6
Noweda	Stemless carline root	12163702	1
Noweda	Tall golden rod, cut	11280306	3
Noweda	Testosterone propionate	1511003	1
Noweda	Testosterone propionate	12163610	3
Noweda	Tetracaine hydrochlor...	120700-02	1
Noweda	Tetracaine hydrochlor...	1106035-01	5
Noweda	Thyme, rubbed	521303908	1
Noweda	Thymol	12k01-no2	1
Noweda	Tonka beans	1040A111026	1
Noweda	Tragacanth, powdered	14162207	1
Noweda	Tretinoin	15F01-B06-315491	2
Noweda	Tretinoin	15078307	1
Noweda	Tretinoin	13D30-N14	3
Noweda	Triamcinolone acetonide	13K15-B31-289029	5
Noweda	Triamcinolone acetonide	12301501	5
Noweda	Triclosan	13247903	2
Noweda	Triclosan	12359107	2
Noweda	Urea	4109M-02544	1
Noweda	Urea	15154110	38
Noweda	Wax, yellow	15023002	1
Noweda	Wax, yellow	150230004	1
Noweda	Wax, yellow	280A120627	7
Noweda	Willow bark, cut	4609A-01305	1
Noweda	Xylometazoline hydroc...	12027902	2
Noweda	Yarrow herb, cut	12378603	1
Noweda	Yarrow herb, cut	13381804	2
Noweda	Zedoary root, cut	11091605	3

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Noweda	Mannitol	1205016-02	8
Noweda / 28.10.13 EuRho	Magnesium sulfate hep...	1211022-01	3
Noweda 06.04.2016	Bitter fennel	15443201	1
Noweda 17.07.2015	Hydrocortisone acetate	14301403	5
Noweda 19.05.2015	Estriol	14272104	2
Noweda/ Fischar	Tartaric acid	702061	1
Noweda/ Infecto Pharm	Salicylic acid	16/15-08-330773	1
Noweda/Anzag	Clobetasol propionate	13269504	8
Noweda/Euro	Neomycin sulfate	1601046,02	1
Noweda/EuRo	Ambroxol hydrochloride	1609002-02	15
Noweda/Fargon	Clotrimazole	17A19-B04-336568	1
Noweda/Fargon	Phenylephrine hydroch...	17C23-B04-336301	4
Noweda/Fargon	Resocinol	16029-B01-329760	1
Noweda/Fargon	Salicylic acid	14L02-B02-303556	1
Noweda/Fargon	Tetracaine base	16/12-B06-329087	1
Noweda/Henry Lamotte	Black mustard seed, p...	1493708	1
Noweda10.02.2017	Tormentil	15364313	1
Noweda21.04.2015	Erythromycin	14044110	2
Noweda29.04.2015	Chloral hydrate	15083701	3
Noweda29.04.2015	Mannitol	14323113	2
Noweda-Fragon	Naphazoline hydrochlo...	1303B03305291	1
NWM CHina	Estriol	ES150304	1
NWM CHina	Pregnenolone	PR150329	1
NWM CHina	Testosterone	TE150402	1
NWM Shandong Northwes...	Prasterone	PST151207	1
NWM Shandong Northwes...	Pregnenolone	PRN151219	1
Oberlinden Apotheke	Betamethasone valerate	c073/115040	1
Oberlinden Apotheke	Betamethasone valerate	C0073/115040	3
Oberlinden Apotheke	Clotrimazole	201401300378	2
Oberlinden Apotheke	Triamcinolone acetonide	B0544/114010	1
Oberlinden- Apotheke	Betamethasone valerate	C0073/1	3
Orgamol Pharma	Tetracaine hydrochlor...	6093547	4
Otto Fischar	Tartaric acid	705061	5
Otto Fischar GmbH&Co...	Citric acid	6494061	11
Otto Fischar GmbH&Co...	Citric acid	6496051	1
Otto Fischer GmbH	Citric acid	6492061	1
otz pharma	Lactose monohydrate	130102802	1
p&m cosmetics/gehe	Clotrimazole	16J07-B02-335773	4
Penny	Sucrose	L12513800	3
Penny Markt	Sucrose	L120612100	3
Penny Markt	Sucrose	201L113530100	3
Penny Markt	Sucrose	L120404100	3
Penny Markt	Sucrose	L120628000	3
Penny Markt	Sucrose	L120325600	2
Peter Kölln	Lactose monohydrate	229719241-0640	1
Pfannenschmidt	Lithium orotate monoh...	ORL00212	1
Pfeifer/Langen	Sucrose	L03144N	1
ph	(Sheng) Di Huang	g211h0532323	1
Ph.Eur.8.0	Glucose monohydrate	1401019-02	1
Phagron	Erythromycin		1
Phagron	Mannitol	1302C004-279405	1
Phagron	Yarrow herb, cut		3
phar	Acetylcysteine	1610011-02	1
Pharmachem	Sucrose	7878	1
Pharmachem	Wax, yellow	16038202	14
Pharmacia & Upjohn Co...	Neomycin sulfate	16165017	1
Pharmacia & Upjohn Co...	Neomycin sulfate	16122623	1
pharm-a-speres	Sucrose	8731003	1
pharm-a-spheres/sana ...	Sucrose	8731306	1
phnoix	Neomycin sulfate	130303301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix	Acetylsalicylic acid	13324107	1
Phoenix	Ambroxol hydrochloride	1604035-02	7
Phoenix	Benzocaine	12H10-N06	5
Phoenix	Betamethasone dipropi...	13236806	3
Phoenix	Budesonide, micronized	14071702	8
Phoenix	Butylated hydroxytolu...	13354609	3
Phoenix	Butylated hydroxytolu...	12382305	1
Phoenix	Butylated hydroxytolu...	3111201	1
Phoenix	Camphor	15316811	7
Phoenix	Cellulose, microcryst...	1309001-02	4
Phoenix	Centaury, cut	11097612	9
Phoenix	Cetostearyl alcohol	14322505	2
Phoenix	Chlorhexidine diacetate	11204408	3
Phoenix	Clindamycin hydrochlo...	1710013-010	1
Phoenix	Clioquinol	16J18-B04-332852	4
Phoenix	Clobetasol propionate	16h22-b07-328012	1
Phoenix	Clobetasol propionate	16H-B07-330169	1
Phoenix	Clobetasol propionate	13K08B30287728	1
Phoenix	Dexamethasone	13g08n13	1
Phoenix	Dexamethasone	13084313	4
Phoenix	Dexamethasone	10369928	1
Phoenix	Dexamethasone	15I04-B04	1
Phoenix	Diltiazem hydrochloride	13346914	6
Phoenix	Erythromycin	15042811	4
Phoenix	Ethacridine lactate m...	100201	2
Phoenix	Fenugreek	11313206	6
Phoenix	Finasteride (form III)	1804010-01	3
Phoenix	Furosemide	17D06-B01-338982	1
Phoenix	Gentamicin sulfate	16b29-b01-328731	1
Phoenix	Gentian root, cut	6311B120104	1
Phoenix	Glucose	15329010	4
Phoenix	Hard fat W35	16k08-b10-332763	1
Phoenix	Hydrocortisone acetate	15A30-B09	1
Phoenix	Hydrocortisone acetate	15J28-B12	2
Phoenix	Hydrocortisone acetate	1568/05131813	1
Phoenix	Hypromellose 4000	13E02-B03	1
Phoenix	Hypromellose 4000	14c05-b03-294708	1
Phoenix	Hypromellose 4000	16B16-B07	1
Phoenix	Iodoform	12070601	3
Phoenix	Juniper		1
Phoenix	Lactose monohydrate	1404PLV-1928	4
Phoenix	Liquorice root, peeled	6381A120611	12
Phoenix	Macrogols 4000	12258901	1
Phoenix	Macrogols 4000	16023905	6
Phoenix	Mannitol	15197314	7
Phoenix	Melatonin	15C20-B01	1
Phoenix	Menthol	16095617	1
Phoenix	Menthol	1403022-02	11
Phoenix	Menthol	16095625	1
Phoenix	Methoxsalen	14F06-B05-296502	1
Phoenix	Metronidazole		2
Phoenix	Metronidazole	16089201	5
Phoenix	Metronidazole	121203002	1
Phoenix	Miconazole nitrate	17C01-B02-334882	2
Phoenix	Miconazole nitrate	15E04-B03-327367	1
Phoenix	Minoxidil	16J25-B05-332092	1
Phoenix	Minoxidil	1602013-02	16
Phoenix	Mometasone furoate, m...	MOM004/0716	1
Phoenix	Oxybutynin hydrochlor...	16/23-B03-328916	1
Phoenix	Oxybutynin hydrochlor...	22032E005	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix	Phenylephrine hydroch. ...	15F09-B08	1
Phoenix	Potassium citrate	1326PLV-1739	1
Phoenix	Potassium citrate	15148708	1
Phoenix	Potassium sorbate	1397235	1
Phoenix	Prednicarbate, micron. ...	16J27-b01-332958	1
Phoenix	Prednisolone, microni. ...	11j26-N04	2
Phoenix	Raspberry leaf, cut	10136409	1
Phoenix	Rhubarb, cut	11057817	3
Phoenix	Ribwort plantain, cut	4271A11016	1
Phoenix	Saccharin sodium	14F30-B02-296322	1
Phoenix	Sage leaf, cut	231C130326	1
Phoenix	Salicylic acid	12H10-No8	1
Phoenix	Sildenafil citrate	1311010-04	6
Phoenix	Sodium cyclamate	787/03100611	1
Phoenix	Sodium edetate	15427704	6
Phoenix	Sodium edetate	16F21-B02-333607	2
Phoenix	Sodium salicylate, po. ...	12158309	1
Phoenix	Sodium tetraborate	1603043-01	1
Phoenix	Sorbic acid	11203813	1
Phoenix	Talcum	1305022-01	1
Phoenix	Tartaric acid	12K29B01	1
Phoenix	Tartaric acid	12044115	1
Phoenix	Tartaric acid	12Ho2-N07	1
Phoenix	Tartaric acid	13116308	3
Phoenix	Tetracaine hydrochlor. ...	15437104	15
Phoenix	Triamcinolone acetoneide	13126-N11	2
Phoenix	Triclosan	17K02-B06-343547	4
Phoenix	Triclosan	12003202	2
Phoenix	Trometamol	16168819	15
Phoenix	Tropicamide	14L04-B02	2
Phoenix	Urea	15154113	24
Phoenix	Wormwood herb, cut	12103003	8
PHOENIX	Quinoline yellow	12042004	1
Phoenix 03.04.2014	Melatonin	12J19-N02	1
Phoenix 04.03.14	Ephedrine hydrochloride	1332303	1
Phoenix 05.07.2014	Guaifenesin	1402023-01	2
Phoenix 05.08.2014	Fructose	FR1336122	3
Phoenix 08.07.2014	Potassium citrate	13317624	1
PHOENIX 10.03.2014	Gentamicin sulfate	1305047	2
PHOENIX 12.05.2014	Sorbic acid	13334008	8
Phoenix 12.12.2016 4. ...	Citric acid	1508808	1
PHOENIX 19.05.2014	Black pepper	13221105	1
Phoenix 20.06.2014	Xylose	13A31-B04-292396	1
Phoenix 22.5.14	Epinephrine hydrogen ...	12337804	1
PHOENIX 24.07.2014	Wax, yellow	13001206	1
PHOENIX 25.04.2014	Chloramphenicol	1308017-02	14
PHOENIX 26.03.2014	Thymol	14025802	5
PHOENIX 29.01.2013	Diphenhydramine hydro. ...	120208	1
Phoenix 29.11.12	Clobetasol propionate	700225	3
PHOENIX 30.05.2014	Glucose	13329514	1
PHOENIX 30.06.2014	Ascorbyl palmitate	UQ21206229	1
PHOENIX 31.01.2014	Tobramycin sulfate	1202006-01	1
PHOENIX am 27.12.2013	Liquorice root, natur. ...	281844	1
Phoenix goettin	Acetylcysteine	1208026-01	1
Phoenix Goettin	Dexamethasone	117113352302	1
Phoenix Goettin	Sodium edetate	13G15-N03	1
Phoenix Goettin	Xylose	13133422	1
Phoenix goettingen	Salicylic acid	14300808	2
Phoenix goettingen 08. ...	Wax, yellow	15098901	2
Phoenix, Jembke	Bismuth gallath, basic	10135101	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phoenix/Cealo	Acetylsalicylic acid	17356301	1
Phoenix/Cealo	Erythromycin	17088004/0918E-01362	1
Phoenix/Cealo	Fusidic acid	701714-0003/0918A-07	1
Phoenix/Cealo	Menthol	1509045-01/4116F-018	1
Phoenix/Cealo	Methacholine chloride	16L20-B01-329968	1
Phoenix/Cealo	Methylcellulose	17170404	1
Phoenix/Cealo	Metronidazole	17058016/0818E-0187	1
Phoenix/Cealo	Salicylic acid	3414M-01045	1
Phoenix/Cealo	Salicylic acid	1733303	1
Phonix	Captopril	14B26-B01	3
Phonix	Ethacridine lactate m...	13I24-N03	8
Phonix	Linseed	13296703	3
Phonix	Naphazoline hydrochlo...	1704005-02	1
Phonix	Poloxamers 407	12/13-N10	1
phönix	Aloes, cape, cut	15283011	1
phönix	Betamethasone valerate	15/16B07-319228	1
phönix	Betamethasone valerate	16e23-b03-327350	1
phönix	Bismuth gallath, basic	15J26-B02-325450	2
phönix	Calcium lactate penta...	15101411	3
phönix	Camphor	444a160817	1
phönix	Chlorhexidine diacetate		1
phönix	Clindamycin hydrochlo...	15392820	1
phönix	Clioquinol	16f18b04328511	1
phönix	Clotrimazole	11305513	8
phönix	Gentamicin sulfate	1214605	1
phönix	Hydrocortisone		2
phönix	Hydroxyethylcellulose...	15e07b05316511	1
phönix	Hydroxyethylcellulose...	1610039-02	9
phönix	Ketamine hydrochloride	14103-B02-311706	1
phönix	Macrogols 4000	1701007-01	3
phönix	Macrogols 4000	1535007	1
phönix	Macrogols 4000	15406302	3
phönix	Macrogols 4000	1605012-01	1
phönix	Macrogols 4000	1502031-02	1
phönix	Magnesium hydrogen ci...	15415015	3
phönix	Magnesium hydrogen ci...	15296710	2
phönix	Methylrosanilinium ch...	13184811	1
phönix	Minoxidil		1
phönix	Mometasone furoate, m...	15B11-b03	1
phönix	Myrrh	3351a160330	1
phönix	Myrrh	15021411	1
phönix	Potassium hydrogen ta...	15354604	3
phönix	Prednisolone acetate	1438910	1
phönix	Progesterone, microni...	16g15-b04-326487	1
phönix	Propranolol hydrochlo...	16e09b09326637	1
phönix	Propranolol hydrochlo...	16e09-b09	1
phönix	Propranolol hydrochlo...	2/AKD3	1
phönix	Quinoline yellow	15276803	3
phönix	Rice starch	84a141114	1
phönix	Sodium alginate	11509036-02	1
phönix	Sodium edetate	13142018	5
phönix	Sucrose	16351701	1
phönix	Tannin	15/AKD65	1
phönix	Tannin	14E07-B03-293744	1
phönix	Tormentil	15364343	1
phönix	Triamcinolone acetone...	15430107	14
phönix	Zinc acetate dihydrate	15383301	1
Phönix	Acetylcysteine	1208026	1
Phönix	Acetylsalicylic acid	13324114	1
Phönix	Arginine	L1302006-01	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Atropine sulfate	15301607	1
Phönix	Bacitracin	c2280065	1
Phönix	Benzocaine	13F06-N10	2
Phönix	Benzoic acid	13074605	4
Phönix	Betamethasone valerate	12229807	1
Phönix	Betamethasone valerate	BV-BLT186/0512	1
Phönix	Bismuth gallath, basic	12-13-N11	1
Phönix	Bitter fennel	12094708	2
Phönix	Calcium gluconate	1606034-01	13
Phönix	Calcium lactate penta...	14H01-B05	1
Phönix	Centaury, cut	11093612	2
Phönix	Chloramphenicol	13A1-802	3
Phönix	Chloramphenicol	13080921	3
Phönix	Chloramphenicol	16i12-b07-327124	1
Phönix	Chloramphenicol	10341710	1
Phönix	Chloramphenicol	10341716	1
Phönix	Chlorhexidine diacetate	Retest09113916	1
Phönix	Chlorhexidine diacetate	9113916	1
Phönix	Chlorhexidine diacetate	9113916	1
Phönix	Cinchocaine hydrochlo...	13E13-N07	1
Phönix	Cinnamon rods (Ceylon...	610A130927	1
Phönix	Citric acid	14186312	18
Phönix	Citric acid monohydrate	1403/3	1
Phönix	Citric acid monohydrate	10111503	1
Phönix	Clindamycin hydrochlo...	15392803	1
Phönix	Clindamycin hydrochlo...	12H30-N04	1
Phönix	Clioquinol	14C21-B08	2
Phönix	Clioquinol	14C21-B06	4
Phönix	Clioquinol	16718b04332852	1
Phönix	Clobetasol propionate	17096706	16
Phönix	Clobetasol propionate		1
Phönix	Cocoa butter	27121	1
Phönix	Coltsfoot leaf, cut	278296	4
Phönix	Common broom herb, cut	12358801	2
Phönix	Creatine	mysjhskeu	1
Phönix	Dandelion herb with r...	6611B130321	29
Phönix	Dequalinium chloride	10127-N02Retest	1
Phönix	Dequalinium chloride	1110002-01	1
Phönix	D-galactose	39223PF3	1
Phönix	Diltiazem hydrochloride	17d11-b11-334837	1
Phönix	Diltiazem hydrochloride	13196106	1
Phönix	Emulsifying cetostear...	1323305	1
Phönix	Ephedrine hydrochloride	13323315	1
Phönix	Equisetum stem, cut	12288807	2
Phönix	Erythromycin	121427016	1
Phönix	Erythromycin	16079205	3
Phönix	Estradiol benzoate	16003707	9
Phönix	Estradiol valerate	15C10-B05-326087	2
Phönix	Gentamicin sulfate	11113309	1
Phönix	Gentian root, cut	6311B120315	8
Phönix	Gentian root, cut	269197	2
Phönix	Ginger, peeled	13310601	1
Phönix	Glucose	13178701	2
Phönix	Glucose monohydrate	1210334-02	1
Phönix	Ground ivy herb, cut	11374608	1
Phönix	Hydrocortisone	11119-N03	1
Phönix	Hydrocortisone	12167010	12
Phönix	Hydrocortisone	13341504	10
Phönix	Hydrocortisone acetate	14144804	6
Phönix	Hydrocortisone acetate	12A13-N02	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Hydrocortisone acetate	12K15-No4	1
Phönix	Hydrocortisone acetate	AKD8	1
Phönix	Hydroxyethylcellulose...	12369711	5
Phönix	Hydroxyethylcellulose...	1304008-02	1
Phönix	Hyssop herb, cut	11374306	1
Phönix	Iceland moss, cut	12187104	3
Phönix	Ivy leaf, cut	11101102	1
Phönix	Jasmine flower	12143206	1
Phönix	Juniper	13020903	1
Phönix	Ketoconazole	12171018	1
Phönix	Ketoconazole	13179602	3
Phönix	Lime flower, cut	2090413	1
Phönix	Lime flower, cut	1611A12120	1
Phönix	Liquorice root, natur...	276425	1
Phönix	Macrogols 1500	14072402	1
Phönix	Macrogols 4000	1707026-01	2
Phönix	Male speedwell wort, ...	13089201	4
Phönix	Mannitol	120501602	2
Phönix	Mannitol	12305723	4
Phönix	Marshmallow root, pee...	2090413	1
Phönix	Marshmallow root, pee...	6071a120314	1
Phönix	Marshmallow root, pee...	6071A130408	14
Phönix	Matricaria flower	1300A120426	1
Phönix	Melatonin	APALC/MCT17030111	1
Phönix	Melatonin	11267111	2
Phönix	Melilot herb, cut	1237903	1
Phönix	Menthol	12246204	2
Phönix	Methyl-4-hydroxybenzo...	2AKD3	1
Phönix	Metronidazole	17058016	15
Phönix	Miconazole nitrate	17217502	16
Phönix	Mometasone furoate, m...	16c24-b03-327346	1
Phönix	Mometasone furoate, m...	13/23-N02	1
Phönix	Motherwort herb, cut	12231503	2
Phönix	Naphazoline hydrochlo...	586550002	1
Phönix	Norethisterone acetate	14H29-B03-315702	1
Phönix	Nystatin	12017005	2
Phönix	Oak bark, cut	10348609	2
Phönix	Peppermint leaf, cut	11327808	2
Phönix	Potassium citrate	1211006-01	2
Phönix	Prednicarbate, micron...	16J27-B01-330976	10
Phönix	Prednisolone acetate	11158008	1
Phönix	Prednisolone hydrogen...	13H24-N04	1
Phönix	Prednisolone, microni...	13A07N14	1
Phönix	Prednisolone, microni...	13E15-N06	17
Phönix	Prednisone	11D08-N05	1
Phönix	Progesterone, microni...	15319604	1
Phönix	Propranolol hydrochlo...	15k11-b01-322717	1
Phönix	Raspberry leaf, cut	12362701	7
Phönix	Rhubarb, cut	12211905	3
Phönix	Rose hips with peel, ...	2841B130523	1
Phönix	Rose hips with peel, ...	10281021	1
Phönix	Rose hips with peel, ...	2841B130523	1
Phönix	Salicylic acid	68/AKD4	1
Phönix	Salicylic acid	11311917	6
Phönix	Salicylic acid	16a18-b04-322905	1
Phönix	Silverweed, cut	11328107	7
Phönix	Sodium benzoate	1510044/01	1
Phönix	Sodium carbonate deca...	13036904	2
Phönix	Sodium edetate	11074702	1
Phönix	Sodium edetate	13E15-B10	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix	Sodium hydrogen carbo...	13186705	4
Phönix	Sodium hydrogen carbo...	13390515	1
Phönix	Sodium monohydrogenph...	1302800,	1
Phönix	Sodium monohydrogenph...	13028101	1
Phönix	Spironolactone	1101013-03	1
Phönix	Stearic acid 50	1224915	1
Phönix	Sulfacetamide	11085910	1
Phönix	Tannin	10313107Retest	1
Phönix	Tartaric acid	11e24-n06	1
Phönix	Testosterone propionate	13K25-B05	2
Phönix	Tetracaine base	3465093	2
Phönix	Tetracaine base	11g12-n13	1
Phönix	Tetracaine base	1303018-01	1
Phönix	Tetracaine hydrochlor...	2015E-02039	1
Phönix	Thymol	11358610	1
Phönix	Tormentil	15364302	1
Phönix	Tretinoin	14D09-B04-293224	2
Phönix	Tretinoin	14D09-B04	2
Phönix	Triamcinolone acetonide	13379902	1
Phönix	Triclosan	13L8-B35	1
Phönix	Triclosan	1516A-03402	1
Phönix	Triclosan	K001.1/14/12/823	1
Phönix	Urea	15087110	24
Phönix	Uridine	721010040	1
Phönix	Wormwood herb, cut	12355806	7
Phönix	Yarrow herb, cut	4181A130513	10
PHÖNIX	Dexamethasone	15K17-B07	1
PHÖNIX	Gentamicin sulfate	15226514	5
Phönix 02.01.15	Melatonin	13H21-N06	1
Phönix 04.03.2014	D-galactose	*6878838-13507	1
Phönix 10.01.2014	Senna leaf, cut	13136707	1
Phönix 14.01.2014	Propyl-4-hydroxybenzo...	11373526	2
phönix 4.12.2013	Methyl-4-hydroxybenzo...	1212011/01	1
Phönix Caesar u. Loretz	Hydrocortisone	13225507	1
Phönix Gotha	Bacitracin	C2280010	1
Phönix Gotha	Betamethasone valerate	14188608	7
Phönix Gotha	Chloral hydrate	13F06-N07	4
Phönix Gotha	Chloral hydrate	14038302	1
Phönix Gotha	Citric acid	13050618	4
Phönix Gotha	Citric acid	13A04-N040049-100	1
Phönix Gotha	Citric acid	12B01-N012125/6000	1
Phönix Gotha	Citric acid	13050609	1
Phönix Gotha	Citric acid	12B01-N011973/6000	1
Phönix Gotha	Citric acid	12B01-N012028/6000	1
Phönix Gotha	Dexamethasone	13352311	2
Phönix Gotha	Griseofulvin	140410	3
Phönix Gotha	Lidocaine	15A15-B01-304284	10
Phönix Gotha	Mometasone furoate, m...	1403015-01	2
Phönix Gotha	Oxytetracycline hydro...	1232814	1
Phönix Gotha	Papaverine hydrochlor...	12166101	1
Phönix Gotha	Papaverine hydrochlor...	13324306	1
Phönix Gotha	Prednisolone, microni...	13343112	4
Phönix Gotha	Procaine hydrochloride	190311	1
Phönix Gotha	Quinolinol sulfate po...	120905	9
Phönix Gotha	Salicylic acid	12374207	6
Phönix Gotha	Sodium acetate trihyd...	12316506	1
Phönix Gotha	Testosterone	1408015-01	1
Phönix Gotha	Trometamol	1319-B02-286158	1
Phönix Gotha	Urea	V096/2012	1
Phönix Gotha	Urea	13064202	13

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phönix WE:04.01.18 E...	Tretinoin	17273102	2
Phönix/Euro	Magnesium chloride he...	1208025-01	2
Phönix/Euro	Polymyxin B sulfate	1505005-010	1
Phönix/Euro	Testosterone propionate	110301002	1
PhönixWE:10.02.2017 ...	Prednisolone acetate	16246806	12
Phyto Comm	Fu Ling	G200H1019522	4
phytocomm	(Ku) Xing Ren	g029h0708322	2
phytocomm	Ce Bai Ye	g050h0975321	2
phytocomm	Chan Tui	g067h1817322	1
phytocomm	Di Long	g148h603321	1
phytocomm	Fu Zi	g005h0801322	1
phytocomm	Lu Jiao Jiao	G323H1158321	1
phytocomm	Niu Bang Zi	g026h0433321	2
phytocomm	Shen Qu	g156h0935421	1
phytocomm	Xian He Cao	g008h0547321	1
phytocomm	Zhi Shi	g044h0922321	1
Phytocomm	(Bai) Dou Kou	G016H0508521	1
Phytocomm	(Bai) Dou Kou	G016H0508521	1
Phytocomm	(Bai) Jiang Can	g269h0513222	1
Phytocomm	(Bai) Jiang Can	g144h1708221	1
Phytocomm	(Bai) Jiang Can	G269H0513321	2
Phytocomm	(Huai) Niu Xi	G003H1905422	2
Phytocomm	(Huai) Niu Xi	G003H1905321	3
Phytocomm	(Huai) Niu Xi	G003H1905222	2
Phytocomm	(Huai) Niu Xi	g003h1905222	1
Phytocomm	(Ku) Xing Ren	g029h0708122	1
Phytocomm	(Ku) Xing Ren	G029H0708125	3
Phytocomm	(Sheng) Di Huang	g211h0532121	1
Phytocomm	(Sheng) Di Huang	G211H0532521	4
Phytocomm	(Sheng) Di Huang	g211h053221	1
Phytocomm	(Shi) Chang Pu	G006H0544122	1
Phytocomm	(Shi) Chang Pu	G006H0544322	1
Phytocomm	(Shi) Chang Pu	g006h0544023	1
Phytocomm	(Shi) Chang Pu	G006HS2700N1	1
Phytocomm	Ai Ye	G033H0613022	1
Phytocomm	Ai Ye	G033H0613421	1
Phytocomm	Bai Bian Dou	g093h0510221	1
Phytocomm	Bai Bian Dou	G093H0510221	1
Phytocomm	Bai Bu	G235H0608321	1
Phytocomm	Bai Guo	G114H0509321	2
Phytocomm	Bai Guo	G114H0509521	1
Phytocomm	Bai He	G139H0607421	3
Phytocomm	Bai He	g139h0607221	1
Phytocomm	Bai Hua She She Cao	G173H1153421	4
Phytocomm	Bai Hua She She Cao	G173H1153222	3
Phytocomm	Bai Jiang Cao	G181H1180422	1
Phytocomm	Bai Jiang Cao	G181H1180321	1
Phytocomm	Bai Mao Gen	G125H0512221	2
Phytocomm	Bai Shao Yao	G179H0502521	5
Phytocomm	Bai Shao Yao	g17h0502323	1
Phytocomm	Bai Shao Yao	g179h0502322	1
Phytocomm	Bai Shao Yao	G179H0502121	1
Phytocomm	Bai Tou Weng	G207H0511521	4
Phytocomm	Bai Xian Pi	G090H0515121	1
Phytocomm	Bai Zhi	G020H0503421	1
Phytocomm	Bai Zhu	G041H0501221	1
Phytocomm	Bai Zhu	G041H0501422	3
Phytocomm	Bai Zhu	g041h0501221	3
Phytocomm	Bai Zi Ren	G052HS187	2
Phytocomm	Bai Zi Ren	g052f110315	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phytocomm	Bai Zi Ren	g952f110315	1
Phytocomm	Ban Lan Gen	G128H0816421	2
Phytocomm	Ban Xia (Jiang)	G191H0526322	1
Phytocomm	Ban Xia (Jiang)	g191h0526121	1
Phytocomm	Ban Zhi Lian	G224H0528221	1
Phytocomm	Ban Zhi Lian	G224H0528521	3
Phytocomm	Bo He	g158h1712022	1
Phytocomm	Bu Gu Zhi	G205H1247421	2
Phytocomm	Bu Gu Zhi	g205h1247222	1
Phytocomm	Bu Gu Zhi	G205H1247222	2
Phytocomm	Cang Er Zi	g251h1401021	1
Phytocomm	Cang Er Zi	G251H1401521	1
Phytocomm	Cang Er Zi	G251H1401322	2
Phytocomm	Cang Zhu	g042h1402221	1
Phytocomm	Ce Bai Ye	G050H0975321	1
Phytocomm	Ce Bai Ye	G050HS245	1
Phytocomm	Chai Hu	G055H1055123	1
Phytocomm	Chai Hu	g055h1055123	1
Phytocomm	Chai Hu	G055H1055422	5
Phytocomm	Chan Tui	G067H1817323	2
Phytocomm	Chan Tui	G067H1817322	3
Phytocomm	Che Qian Zi	410704805	1
Phytocomm	Che Qian Zi	G192H0723421	2
Phytocomm	Chen Pi	g073h1128222	1
Phytocomm	Chen Pi	G073H1128222	2
Phytocomm	Chen Pi	g07h1128221	1
Phytocomm	Chi Shao (Yao)	G180H0705422	1
Phytocomm	Chi Shao (Yao)	G189H0705321	1
Phytocomm	Chi Shao (Yao)	G180H0705321	1
Phytocomm	Chi Shao (Yao)	g180h0705122	1
Phytocomm	Chong Wei Zi	g317bp021046	1
Phytocomm	Chuan Bei Mu	G107BP020722A	1
Phytocomm	Chuan Bei Mu	g107f110721	2
Phytocomm	Chuan Lian Zi	G157H0306421	3
Phytocomm	Chuan Mu Tong	G009H0310321	1
Phytocomm	Chuan Mu Tong	g009h0310321	1
Phytocomm	Chuan Niu Xi	g318h0431321	1
Phytocomm	Chuan Xiong	G136H0305421	2
Phytocomm	Chuan Xiong	g136h0305121	1
Phytocomm	Chuang Mu Xiang	g219h0309322	1
Phytocomm	Chuang Mu Xiang	G219H0309322	1
Phytocomm	Chuang Mu Xiang	g219h0309321	1
Phytocomm	Chuang Mu Xiang	G219H0309521	2
Phytocomm	Da Fu Pi	G293H0321221	1
Phytocomm	Da Huang	G212H0320521	3
Phytocomm	Da Huang	g212h0320023	1
Phytocomm	Da Huang	G212H0320321	3
Phytocomm	Da Qing Ye	G127H0318421	2
Phytocomm	Da Zao	g129h0316221	1
Phytocomm	Da Zao	G129H0316521	3
Phytocomm	Dan Dou Chi	G232HS239	1
Phytocomm	Dan Shen	g214h0441121	1
Phytocomm	Dan Shen	G214H0441121	1
Phytocomm	Dan Shen	G214H0441322	2
Phytocomm	Dan Shen	G214H0441522	2
Phytocomm	Dan Zhu Ye	G046H1105321	2
Phytocomm	Dang Gui	g022h1306121	3
Phytocomm	Dang Gui	G022H1306422	3
Phytocomm	Dang Shen	G077H2001521	4
Phytocomm	Dang Shen	G077H2001022	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phytocomm	Deng Xin Cao	f100707	1
Phytocomm	Di Fu Zi	G131H0605122	1
Phytocomm	Di Gu Pi	g149h060112	1
Phytocomm	Di Gu Pi	G149H0601421	3
Phytocomm	Di Long	G148H0603421	2
Phytocomm	Di Long	G148H0603321	2
Phytocomm	Di Yu	G215H0604322	2
Phytocomm	Ding Xiang	G058H0210021	1
Phytocomm	Dong Gua Zi	G049H0551221	1
Phytocomm	Du Huo	g021h610221	1
Phytocomm	Du Huo	G021H1610421	2
Phytocomm	Du Zhong	G101H0736321	1
Phytocomm	Du Zhong	g101h0736021	1
Phytocomm	Du Zhong	G101H0736522	5
Phytocomm	Du Zhong	g101h076521	1
Phytocomm	E Jiao	g036h0861221	1
Phytocomm	E Jiao	G036H0861421	2
Phytocomm	E Jiao	g036h0861121	1
Phytocomm	E Zhu	h1110021	1
Phytocomm	E Zhu	H1110021	1
Phytocomm	Fang Feng	g217h0728923	1
Phytocomm	Fang Feng	g217h0728221	1
Phytocomm	Fang Feng	G217H0728521	5
Phytocomm	Fang Ji	G236H0727321	1
Phytocomm	Fu Ling	G200H1019321	2
Phytocomm	Fu Ling	g200h1019321	3
Phytocomm	Fu Pen Zi	g213h1813122	1
Phytocomm	Fu Shen	G201H1021422	2
Phytocomm	Fu Shen	G201H1021321	1
Phytocomm	Fu Shen	g201h1021321	1
Phytocomm	Fu Zi	G005H0801322	1
Phytocomm	Fu Zi	g005h0801122	1
Phytocomm	Fu Zi	G005H0801524	1
Phytocomm	Gan Cao	g119h0521121	1
Phytocomm	Gan Cao	G119H0521121	5
Phytocomm	Gan Cao	G119H0521522	8
Phytocomm	Gan Jiang	G252H1145221	1
Phytocomm	Gan Jiang	G252H1145421	1
Phytocomm	Gan Jiang	g252h01145221	1
Phytocomm	Gan Jiang	G252H1145521	1
Phytocomm	Ge Gen	G206H1301221	2
Phytocomm	Ge Gen	G206H1301422	2
Phytocomm	Gou Qi Zi	G150H0913422	4
Phytocomm	Gou Qi Zi	g150h0913221	1
Phytocomm	Gou Teng	G247H1210422	2
Phytocomm	Gou Teng	g247h1210122	1
Phytocomm	Gu Sui Bu	g094h1011023	1
Phytocomm	Gu Sui Bu	g094h1011221	1
Phytocomm	Gu Sui Bu	G094H1011221	3
Phytocomm	Gua Lou	G241H1096521	1
Phytocomm	Gua Lou	G241H1096321	1
Phytocomm	Gua Lou	g241h1096321	1
Phytocomm	Gua Lou Ren	g243h1095123	1
Phytocomm	Guang Huo Xiang	G007H2005521	2
Phytocomm	Gui Zhi	g070h1001422	1
Phytocomm	Gui Zhi	G070H1001221	2
Phytocomm	Gui Zhi	G070H1001523	2
Phytocomm	Gui Zhi	G070H100422	1
Phytocomm	Han Lian Cao	G095H0770321	3
Phytocomm	Han Lian Cao	g095h0770321	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phytocomm	Han Lian Cao	G095H0770521	4
Phytocomm	Hard fat in pastilles	209314	1
Phytocomm	He Huan Pi	G010H0635423	1
Phytocomm	He Huan Pi	g010h0635221	1
Phytocomm	He Shou Wu	G198H0718222	2
Phytocomm	He Shou Wu	G198H0718421	4
Phytocomm	He Shou Wu	g198h0718222	1
Phytocomm	Hong Hua	G057H0901221	1
Phytocomm	Hong Hua	g057h0901121	1
Phytocomm	Hong Jing Tian	G237H0906321	1
Phytocomm	Hong Jing Tian	G237H0906221	1
Phytocomm	Hou Po	G154H0917222	1
Phytocomm	Hu Zhang	G197H0851321	1
Phytocomm	Hu Zhang	g197h0851321	1
Phytocomm	Hua Shi	g238h1332121	1
Phytocomm	Huai Hua	G263H1423421	1
Phytocomm	Huang Bo	h1203022	1
Phytocomm	Huang Bo	g188h1203123	1
Phytocomm	Huang Jing	g196h1206121	1
Phytocomm	Huang Jing	G196H1206121	1
Phytocomm	Huang Lian	g079h1205321	1
Phytocomm	Huang Lian	G079H1205421	1
Phytocomm	Huang Lian	G079H1205521	4
Phytocomm	Huang Lian	G079H1205321	6
Phytocomm	Huang Qi	g040h1202321	2
Phytocomm	Huang Qi	G040H1202423	2
Phytocomm	Huang Qi	G040H1202521	2
Phytocomm	Huang Qin	G225H1204522	3
Phytocomm	Huang Qin	g225h1204121	1
Phytocomm	Huang Qin	G225H1204221	2
Phytocomm	Huo Ma Ren	G141HS061	2
Phytocomm	Ji Li	G240H1410421	1
Phytocomm	Ji Nei Jin	g259h2102321	1
Phytocomm	Ji Nei Jin	G259H2102321	1
Phytocomm	Ji Xue Teng	g159h2101122	1
Phytocomm	Ji Xue Teng	G159H2101423	3
Phytocomm	Jiang Huang	G085H1716422	4
Phytocomm	Jiao Gu Lan	G065FB1276421	3
Phytocomm	Jie Cao	H2111921	1
Phytocomm	Jie Geng	g193h1015221	1
Phytocomm	Jie Geng	G193H1015521	4
Phytocomm	Jin Qian Cao	g151h0847122	1
Phytocomm	Jing Jie	g221h1038221	1
Phytocomm	Jing Jie	G221H1038221	1
Phytocomm	Ju Hua	G064H1231421	3
Phytocomm	Jue Ming Zi	G059H1097221	3
Phytocomm	Ku Shen	G233H0908121	1
Phytocomm	Ku Shen	h0908121	1
Phytocomm	Lai Fu Zi	G209H1220221	1
Phytocomm	Lai Fu Zi	G209H1220321	1
Phytocomm	Lian Qiao	G105H1132221	1
Phytocomm	Lian Qiao	G105H1132521	1
Phytocomm	Ling Zhi	g109fb2405321	1
Phytocomm	Long Dan (Cao)	G112H1705321	3
Phytocomm	Long Dan (Cao)	g112h1705022	1
Phytocomm	Lu Gen	G189H2014221	1
Phytocomm	Lu Gen	G189H2014521	3
Phytocomm	Lu Jiao Shuang	FB1419901	1
Phytocomm	Lu Jiao Shuang	G062H1419421	1
Phytocomm	Lu Lu Tong	G142H1336521	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phytocomm	Lu Lu Tong	G142H1336321	1
Phytocomm	Ma Bian Cao	G248H1085221	1
Phytocomm	Ma Huang	G097H1101522	2
Phytocomm	Mai Men Dong	G176H1118521	3
Phytocomm	Mai Men Dong	g176h1118221	1
Phytocomm	Mai Ya	g123h1119022	1
Phytocomm	Mai Ya	G123H1119521	1
Phytocomm	Mi Huan Jun	G510H0418422	1
Phytocomm	Mi Huan Jun	G510H0418221	1
Phytocomm	Mo Yao	G1690732222	1
Phytocomm	Mu Dan Pi	g166j074122	1
Phytocomm	Mu Dan Pi	G166H0741321	3
Phytocomm	Mu Dan Pi	G166H0741122	2
Phytocomm	Mu Dan Pi	G166H0741121	1
Phytocomm	Mu Dan Pi	G166H0741421	4
Phytocomm	Mu Gua	G063H0424521	4
Phytocomm	Mu Zei	G265H0427421	1
Phytocomm	Niu Bang Zi	G026H0433521	3
Phytocomm	Niu Bang Zi	G026H0433321	3
Phytocomm	Nü Zhen Zi	420319901	1
Phytocomm	Nü Zhen Zi	G138H0345421	2
Phytocomm	Pi Pa Ye	G099H0832221	1
Phytocomm	Pu Gong Ying	G239H1315321	6
Phytocomm	Qian Hu	G185H0926322	4
Phytocomm	Qian Shi	G102H0836422	1
Phytocomm	Qiang Huo	G172H0840321	2
Phytocomm	Qiang Huo	g172h0841321	1
Phytocomm	Qing Pi	g074h100713	1
Phytocomm	Quan Xie	G222F110314	1
Phytocomm	Ren Shen	g115ho562121	1
Phytocomm	Rou Gui	g069h0617123	1
Phytocomm	Rou Gui	H0617022	1
Phytocomm	Rou Gui	G069H0617323	2
Phytocomm	Ru Xiang	G174H0820421	3
Phytocomm	Ru Xiang	G174H0820231	1
Phytocomm	San Leng	g270h0311221	1
Phytocomm	San Qi	g203h0312021	1
Phytocomm	Sang Bai Pi	G161H1042321	3
Phytocomm	Sang Ji Shend	G146H1045522	4
Phytocomm	Sang Ji Shend	g146h1045222	1
Phytocomm	Sang Piao Xiao	G282H421029201	1
Phytocomm	Sang Shen	G163H1047121	1
Phytocomm	Sang Shen	G163H1047221	2
Phytocomm	Sang Ye	G162H1043121	1
Phytocomm	Sang Ye	G162H1043421	4
Phytocomm	Sha Ren	g018h0939321	1
Phytocomm	Sha Ren	G018H0939522	1
Phytocomm	Sha Ren	G018H0939321	1
Phytocomm	Sha Shen (Bei)	G117H0754322	4
Phytocomm	Shan Yao	G091H0330521	4
Phytocomm	Shan Yu Rou	g080h0333221	1
Phytocomm	Shan Yu Rou	G080H0333221	1
Phytocomm	Shan Yu Rou	G080H0333421	1
Phytocomm	Shan Zha	G082H0331521	1
Phytocomm	Shan Zha	G082H0331324	1
Phytocomm	She Chuang Zi	G076H1154122	1
Phytocomm	She Gan	G048H1034422	1
Phytocomm	Shen Qu	G156H0935521	2
Phytocomm	Sheng Jiang	G253H0534221	1
Phytocomm	Sheng Jiang	G253H0534522	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phytocomm	Sheng Ma	g068h04317121	1
Phytocomm	Sheng Ma	G068H0437421	2
Phytocomm	Shi Gao	G120H0541223	3
Phytocomm	Shu Di (Huang)	G210H1601521	3
Phytocomm	Shu Di (Huang)	G210H1601421	6
Phytocomm	Shu Di (Huang)	G210H1601221	1
Phytocomm	Suan Zao Ren	G254H1418421	1
Phytocomm	Suan Zao Ren	g254h1418121	2
Phytocomm	Tai Zi Shen	G204H0446521	1
Phytocomm	Tai Zi Shen	G204H0446221	1
Phytocomm	Tao Ren	g184h1007122	1
Phytocomm	Tian Hua Fen	h0417022	1
Phytocomm	Tian Hua Fen	G242H0417423	2
Phytocomm	Tian Hua Fen	G242H0417522	2
Phytocomm	Tian Men Dong	G037H0415321	1
Phytocomm	Tian Men Dong	g03h0415521	1
Phytocomm	Tian Men Dong	G037H0415221	1
Phytocomm	Tian Men Dong	g037h0415221	1
Phytocomm	Tu Fu Ling	G231H0325421	3
Phytocomm	Tu Fu Ling	G231H0325121	2
Phytocomm	Tu Si Zi	G086H1243521	4
Phytocomm	Tu Si Zi	h1243022	1
Phytocomm	Wang Bu Liu Xing	G281H0456521	1
Phytocomm	Wei Ling Xian	G075H0943422	1
Phytocomm	Wei Ling Xian	g075h0943223	1
Phytocomm	Wei Ling Xian	G075H0943223	3
Phytocomm	Wu Gong	G311H1411321	1
Phytocomm	Wu Gong	G311H411321	1
Phytocomm	Wu Jia Pi	g002h0402321	1
Phytocomm	Wu Mei	G167H1026421	1
Phytocomm	Wu Mei	G167H1026221	1
Phytocomm	Wu Yao	g140h1025221	1
Phytocomm	Wu Yao	G140H1025521	3
Phytocomm	Wu Zhu Yu	G103H0766022	1
Phytocomm	Xi Xian Cao	G229H1431321	2
Phytocomm	Xi Yang Shen	G324HS4070L1	1
Phytocomm	Xi Yang Shen	g324h0647321	1
Phytocomm	Xi Yang Shen	G324H0647321	1
Phytocomm	Xia Ku Cao	G202H1064221	4
Phytocomm	Xia Ku Cao	g202h1064121	1
Phytocomm	Xia Ku Cao	G202H1064421	1
Phytocomm	Xian He Cao	G008H0547321	1
Phytocomm	Xian He Cao	G008H0547121	1
Phytocomm	Xian Mao	G304H0549321	2
Phytocomm	Xian Mao	G304H0549521	2
Phytocomm	Xiang Fu	g088h0930123	1
Phytocomm	Xiang Fu	G088H0930521	6
Phytocomm	Xiao Hui Xiang	G104H0339522	3
Phytocomm	Xiao Hui Xiang	g104h0339321	1
Phytocomm	Xie Bai	G504H1830421	1
Phytocomm	Xie Bai	G504H1830422	1
Phytocomm	Xin Yi	G153H0758521	1
Phytocomm	Xu Duan	G092H2106422	1
Phytocomm	Xuan Fu Hua	G126H1168521	1
Phytocomm	Xuan Shen	g22h0537322	1
Phytocomm	Xuan Shen	G223H0537421	1
Phytocomm	Yan Hu Suo	g081f1000301	1
Phytocomm	Yan Hu Suo	G081H0805521	2
Phytocomm	Ye Jiao Teng	G288H0885223	1
Phytocomm	Ye Jiao Teng	g228h0085223	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Phytocomm	Ye Jiao Teng	G288H0885522	4
Phytocomm	Yi Mu Cao	G134H1050521	4
Phytocomm	Yi Mu Cao	g134h1050024	1
Phytocomm	Yi Yi Ren	G078H1701522	3
Phytocomm	Yi Yi Ren	g078h1701323	1
Phytocomm	Yi Yi Ren	G078H1701323	7
Phytocomm	Yi Zhi Ren	G015H1051521	1
Phytocomm	Yi Zhi Ren	G015H1051321	2
Phytocomm	Yin Yang Huo	G098H1114422	2
Phytocomm	Yu Xing Cao	g124h1149122	1
Phytocomm	Yu Xing Cao	G124H1149423	2
Phytocomm	Yu Xing Cao	g124h1149222	1
Phytocomm	Yu Xing Cao	G124H1149222	1
Phytocomm	Yu Zhu	G195H0517421	2
Phytocomm	Yuan Zhi	g194h140622	1
Phytocomm	Yuan Zhi	G194H1406222	2
Phytocomm	Yuan Zhi	G194H1406522	3
Phytocomm	Ze Lan	g268h1606221	1
Phytocomm	Ze Lan	G268H1606421	1
Phytocomm	Ze Xie	g011h1605122	1
Phytocomm	Ze Xie	G011H1605422	1
Phytocomm	Ze Xie	G011H1605521	6
Phytocomm	Zhe Bei Mu	G108H0712422	1
Phytocomm	Zhen Zhu Mu	g155h0979221	1
Phytocomm	Zhi Gan Cao	g118h0881322	1
Phytocomm	Zhi Gan Cao	G118H0881522	4
Phytocomm	Zhi Ke	G043H0921422	2
Phytocomm	Zhi Ke	G043H0921221	1
Phytocomm	Zhi Shi	G044H0922521	3
Phytocomm	Zhi Shi	g044h0922221	1
Phytocomm	Zhi Shi	G044H0922221	5
Phytocomm	Zhi Zi	G110H1137521	4
Phytocomm	Zhi Zi	g110h1137121	1
Phytocomm	Zhu Ling	G199H1509421	2
Phytocomm	Zhu Ru	G047H0625321	3
Phytocomm	Zhu Ru	g047h0625321	2
Phytocomm	Zhu Ru	G047h625321	1
Phytocomm	Zhu Ru	G047H0625422	3
Phytocomm	Zi Wan	G038H1215321	2
PhytoComm	(Bai) Jiang Can	G269H0513222	1
PhytoComm	(Fen) Bi Xie	G0290H1251121	1
PhytoComm	(Huai) Niu Xi	G003H190522	1
PhytoComm	(Ku) Xing Ren	G029H0708322	2
PhytoComm	(Ku) Xing Ren	G029H0708622	1
PhytoComm	(Ku) Xing Ren	G029H708125	1
PhytoComm	(Sheng) Di Huang	G211H0532621	1
PhytoComm	Ai Ye	G033H0613021	1
PhytoComm	Bai Bu	G235H0608221	1
PhytoComm	Bai He	G139H0607221	1
PhytoComm	Bai Shao Yao	G179H0502322	2
PhytoComm	Bai Shao Yao	1214401	1
PhytoComm	Bai Tou Weng	G207H0511221	3
PhytoComm	Bai Zhu	G041H0501621	1
PhytoComm	Bie Jia	G507H2401321	1
PhytoComm	Bie Jia	G507H2401121	1
PhytoComm	Bo He	G158H1712121	1
PhytoComm	Cao Guo	E241050	1
PhytoComm	Chai Hu	G055H1055322	1
PhytoComm	Che Qian Zi	G192H0723321	2
PhytoComm	Che Qian Zi	G192H0723621	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PhytoComm	Chuan Bei Mu	G107F110721	1
PhytoComm	Chuan Mu Tong	G009H0310521	1
PhytoComm	Chuan Niu Xi	G318H0431321	1
PhytoComm	Chuan Xiong	G136H0305221	3
PhytoComm	Dan Shen	G214H0441221	3
PhytoComm	Dang Gui	G022H1306621	4
PhytoComm	Dang Gui	G022H1306121	1
PhytoComm	Dang Gui	G13-0784	1
PhytoComm	Dang Shen	G077H2001621	1
PhytoComm	Di Long	G148H0603022	1
PhytoComm	E Bu Shi Cao	G061H180122	1
PhytoComm	E Jiao	G036H0861221	4
PhytoComm	E Zhu	G083H1110222	2
PhytoComm	Fang Ji	G236H0727123	1
PhytoComm	Fu Ling	G200H1019622	1
PhytoComm	Fu Ling	G200H1019221	2
PhytoComm	Gan Cao	G119H0521524	1
PhytoComm	Ge Gen	G206H1301323	1
PhytoComm	Geng Mi	G329H1320121	1
PhytoComm	Gou Qi Zi	G150H0913221	1
PhytoComm	Gu Sui Bu	G094H1011522	1
PhytoComm	Gua Lou Ren	G243H1095421	1
PhytoComm	He Zi	G291H1239221	2
PhytoComm	Huang Bo	G188H1203321	4
PhytoComm	Huang Qi	G040H1202621	1
PhytoComm	Huang Qi	G040H1202321	2
PhytoComm	Huang Qin	G225H1204322	5
PhytoComm	Jiang Huang	G085H1716321	1
PhytoComm	Ju Hua	G064H1231521	1
PhytoComm	Jue Ming Zi	G059H1097221	1
PhytoComm	Kuan Dong Hua	H1219921	1
PhytoComm	Lian Qiao	G105H1132321	1
PhytoComm	Ling Zhi	G109F-B2405221	1
PhytoComm	Long Gu	G025H1706122	1
PhytoComm	Long Yan Rou	G144H1708221	1
PhytoComm	Lu Lu Tong	G142H1336122	1
PhytoComm	Ma Chi Xian	G255F110310	1
PhytoComm	Man Jing Zi	G250H1505122	1
PhytoComm	Mo Yao	G169H0732222	2
PhytoComm	Mu Gua	G063H0424221	2
PhytoComm	Mu Zei	G265H0427823	1
PhytoComm	Qian Shi	G102H0836121	1
PhytoComm	Qin Jiao	G113H1059121	1
PhytoComm	Qin Jiao	G113H1059521	1
PhytoComm	Qing Pi	G074H0809621	1
PhytoComm	Rou Gui	G069H0617123	2
PhytoComm	San Leng	G270H0311221	1
PhytoComm	Sang Bai Pi	G161H1042322	1
PhytoComm	Sang Ji Shend	G146H1045322	1
PhytoComm	Sha Ren	G018H0939523	1
PhytoComm	Shan Yu Rou	G080H0333621	1
PhytoComm	She Chuang Zi	G076H1154021	1
PhytoComm	Shen Qu	G156H0935522	1
PhytoComm	Shen Qu	G156H0935121	1
PhytoComm	Sheng Jiang	G253H0534321	1
PhytoComm	Shu Di (Huang)	G210H1601121	1
PhytoComm	Suo Yang	G260H1831121	1
PhytoComm	Tao Ren	G184H1007321	1
PhytoComm	Tu Fu Ling	G23H0325121	1
PhytoComm	Wang Bu Liu Xing	G281H0456322	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
PhytoComm	Wang Bu Liu Xing	G281H456421	1
PhytoComm	Wang Bu Liu Xing	G281H0456221	1
PhytoComm	Wei Ling Xian	G075H0943223	1
PhytoComm	Wu Jia Pi	G002H0402321	4
PhytoComm	Wu Yao	G140H1025221	2
PhytoComm	Wu Zhu Yu	G103H0766321	1
PhytoComm	Xi Xian Cao	G229H1431521	1
PhytoComm	Xiang Fu	G088H0930321	1
PhytoComm	Xie Bai	G504H1830321	1
PhytoComm	Xin Yi	G153H0758321	1
PhytoComm	Yan Hu Suo	G081H0805223	2
PhytoComm	Yan Hu Suo	G081H0805223	1
PhytoComm	Yi Mu Cao	G134H1050321	1
PhytoComm	Yi Mu Cao	G134H1050621	1
PhytoComm	Yu Zhu	G195H0517121	1
PhytoComm	Ze Lan	G268H1606221	2
PhytoComm	Ze Xie	G011H1605122	1
PhytoComm	Zhe Bei Mu	G108H0712321	1
PhytoComm	Zhi Gan Cao	G118H0881322	1
PhytoComm	Zhi Mu	G019H0824321	1
PhytoComm	Zhi Shi	G044H0922321	4
PhytoComm	Zhi Zi	G110H1137121	1
PhytoComm	Zhu Ling	G199H1509521	1
PhytoComm	Zhu Ling	G199H1509022	1
PhytoComm	Zi Hua Di Ding	G249H1216321	2
PhytoComm	Zi Su Ye	G182H1214321	3
PhytoComm	Zi Su Ye	G182HS2680P1	1
Polpharma/Phoenix	Sildenafil citrate	131101001	1
Post-Apotheke 22.01.2...	Camphor	11032012	1
Probe	Elder root, cut	12345	1
Probe	Peppermint leaf, cut	1223456	1
Purren Apotheke	Carnosine (L-) (H1014)	1813979202	1
Purren Apotheke	Chaga (H1011)	1514675601	1
Purren Apotheke	Chondroitin sulfate ...	1821964101	2
Purren Apotheke	Creatine (H1031)	965201	1
Purren Apotheke	Glutathione (H1020)	1254955501	1
Purren Apotheke	Hard fat in pastilles	10I03-N08	1
Purren Apotheke	Heavy alkaline bismut...	11H31-N06	1
Purren Apotheke	Inositol hexanicotina...	1811979201	2
Purren Apotheke	Isoflavones (H1025)	1314955701	1
Purren Apotheke	Leucine (L-) (H1046)	1809965101	2
Purren Apotheke	Ornithine (L-) (H1048)	1817979301	1
Purren Apotheke	Threonine (L-) (H1052)	1833014901	1
Purren Apotheke	Tryptophan (L-) (H1053)	1825950001	1
Purren Apotheke	Tryptophan (L-) (H1053)	978501	1
Purren Apotheke	Valine (L-) (H1056)	1807936701	1
Purren Apotheke	Vitamin B6 (H1103)	1554938001	1
Purren Apotheke	Vitamin E NAT (H1107)	945901	1
Raben-Apotheke Mildenau	Clobetasol propionate	13k08-b30-287727	1
Raben-Apotheke Mildenau	Dequalinium chloride	1405012-01	4
Raben-Apotheke Mildenau	Urea	13J17-Bo2-287604	1
Reagenz	Caffeine		1
Redwood	Damiana leaf, cut	216451	1
Rezepturlager - Verlä...	Carbomers 35000 / 50000	115/14	1
Rezepturlager - Verlä...	Ethacridine lactate m...	121/15	1
Rezepturlager - Verlä...	Ethacridine lactate m...	220/17	1
Rezepturlager - Verlä...	Hydroxyethylcellulose...	54/16	1
Rhodia GmbH/Noweda	Vanillin	12D02-N03	1
Richter	Common wormwood herb,...	b200112	1
Richter	Elecampane rhizome, cut	P158683	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Richter	Gentian root, cut	P200907	1
Richter	Ginger, peeled	KR27992	1
Richter	Ground ivy herb, cut	s99/09	1
Richter	Horse chestnut flower...	s-09/09	1
Richter	Lungwort herb, cut	P195884	1
Richter	Maize style, cut	B189316	1
Richter	Red-rose petal, cut	P236157	1
Richter	Tormentil	P209904	1
Rondell Apotheke	Caffeine citrate	15A05-B03-307309	2
Rondell Apotheke	Glucose	15329003	1
Rondell Apotheke	Lidocaine	15332304	1
Rondell Apotheke	Polymyxin B sulfate	121161	1
Rondell Apotheke	Tetracaine base	1312034-01	2
Rondell Apotheke	Triclosan	11965460	1
s. Protokoll	Tetracycline hydrochl...	G270416A(int.Ch.Be	1
S.Binder	Dequalinium chloride	1210008-04	1
s.Prüfprotokoll	Benzocaine	s.Prüfprotokoll	1
s.Prüfprotokoll	Hydroxyethylcellulose...	s.Prüfprotokoll	1
s.Prüfprotokoll	Hydroxyethylcellulose...	s.Prüfprotokoll	1
s.Prüfprotokoll	Potassium sorbate	s.Prüfprotokoll	1
SA	Clotrimazole	14C04-B30-291661	3
SA	Oxytetracycline hydro...	140701201	1
Sana 04.02.17 Ek:13,83	Salicylic acid	16010910	15
sana 09.12.14 3,75	Guaiacum wood, cut	13343802	1
Sana 21.10.14 Ek 3.57	Maize style, cut	13381603	2
Sanacorp	Magnesium hydrogen ci...	13389108	2
Sanacorp	Diltiazem hydrochloride	14059523	1
Sanacorp	Histamine dihydrochlo...	15C25-B01-306047	1
Sanacorb	Potassium citrate	12135311	3
Sanacorb;20.04.15; 7,...	Citric acid	1375410	1
sanacorp	Bismuth gallath, basic	13/1818B13-292057	1
sanacorp	Calcium gluconate	1407041-01	7
sanacorp	Caraway fruit	2656259	1
sanacorp	Melatonin	14110320	1
sanacorp	Oxalic acid	11098213	1
sanacorp	Prednisolone, microni...	14/10-B03-302266	3
sanacorp	Prednisone	13262302	1
sanacorp	Progesterone, microni...	15124902	2
sanacorp	Pumpkin seed, peeled	14186909	1
sanacorp	Quinoline yellow	3214M-03146	1
sanacorp	Rice starch	94A130304	2
sanacorp	Sodium edetate	16F21-B02-324370	4
Sanacorp	Agnus castus fruit	12274006	3
Sanacorp	Alchemilla, cut	12153309	14
Sanacorp	Alfatradiol	13201605	4
Sanacorp	Allantoin	12374409	19
Sanacorp	Aluminium potassium s...	12184402	2
Sanacorp	Ambroxol hydrochloride	18013219	3
Sanacorp	Ambroxol hydrochloride	170204-02	1
Sanacorp	Angelica root	277310	1
Sanacorp	Bearberry leaf, cut	2391B110209	3
Sanacorp	Benzocaine	13213506	1
Sanacorp	Betamethasone dipropi...	17105911	5
Sanacorp	Betamethasone dipropi...	BDP102/0817	1
Sanacorp	Betamethasone dipropi...	1705020-01	18
Sanacorp	Betamethasone valerate	CH.-B.AKD73	1
Sanacorp	Betony herb, cut	13048903	3
Sanacorp	Bilberry leaf, cut	2191A130410	2
Sanacorp	Blackcurrant leaf, cut	2251A110204	4
Sanacorp	Boric acid	1/AKD3	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Calendula flower, cut	12221601	4
Sanacorp	Camphor	14310720	2
Sanacorp	Camphor	17004306	9
Sanacorp	Camphor	18019802	6
Sanacorp	Camphor	444A130115	1
Sanacorp	Cannabidiol	H180.119-2	16
Sanacorp	Cannabidiol	CBDAPI1804	11
Sanacorp	Caraway fruit	2770A120905	4
Sanacorp	Caraway fruit		2
Sanacorp	Carbomers 35000 / 50000	12274922	2
Sanacorp	Carbomers 35000 / 50000	6/AKD13	1
Sanacorp	Carbomers 35000 / 50000	160702801	1
Sanacorp	Carmellose sodium (40...	AKD6	1
Sanacorp	Cetostearyl alcohol	575778	1
Sanacorp	Cetyl palmitate 15	8/AKD13	1
Sanacorp	Cetyl palmitate 15	8/AKD12	1
Sanacorp	Chicory herb, cut	12369101	6
Sanacorp	Chloral hydrate	17346301	12
Sanacorp	Chloramphenicol	41/AKD23	1
Sanacorp	Chloramphenicol	41/AKD24	1
Sanacorp	Chlorhexidine diacetate	18072403	1
Sanacorp	Chlorhexidine diacetate	7/AKD6	1
Sanacorp	Chlorhexidine diacetate	CH.B.6AKD/1	1
Sanacorp	Ciclopirox olamine	17118806	11
Sanacorp	Cinnamon rods (Ceylon...	13231601	1
Sanacorp	Citric acid	12156811	11
Sanacorp	Citric acid	16208608	5
Sanacorp	Citric acid	40/AKD2	1
Sanacorp	Citric acid	12156832	8
Sanacorp	Citric acid	1216832	1
Sanacorp	Citric acid	41/AKD4	1
Sanacorp	Citric acid monohydrate	1215937	1
Sanacorp	Clindamycin hydrochlo...	K-171219B	1
Sanacorp	Clindamycin hydrochlo...	16D05-B08-332106	2
Sanacorp	Clindamycin hydrochlo...	18D11-B01-352732	1
Sanacorp	Clindamycin hydrochlo...	1511004-03	3
Sanacorp	Clindamycin hydrochlo...	16H03-B02-341602	1
Sanacorp	Clindamycin hydrochlo...	16H03-B02-339275	2
Sanacorp	Clioquinol	11304111	4
Sanacorp	Clioquinol	8/AKD4	1
Sanacorp	Clobetasol propionate	12107608	4
Sanacorp	Clobetasol propionate	13g23n01	1
Sanacorp	Clotrimazole	18147212	4
Sanacorp	Clotrimazole	1230805	1
Sanacorp	Clotrimazole	17360203	7
Sanacorp	Clotrimazole	15g14-b05-311100	1
Sanacorp	Clotrimazole	105/AKD1	1
Sanacorp	Clotrimazole	13224214	10
Sanacorp	Clotrimazole	13184105	15
Sanacorp	Clotrimazole	18147102	5
Sanacorp	Cocoa butter	AKD27	1
Sanacorp	Cocoa butter	AKD28	1
Sanacorp	Cocoa butter	17055908	8
Sanacorp	Cocoa butter	170559908	1
Sanacorp	Cyproterone acetate	18B05-B05-346307	1
Sanacorp	Dandelion herb with r...	13108803	2
Sanacorp	Diclofenac sodium	3616Q-02994	1
Sanacorp	Diltiazem hydrochloride	11134121	2
Sanacorp	Diltiazem hydrochloride	11134115	4
Sanacorp	Diltiazem hydrochloride	14c12-b04-292757	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Diltiazem hydrochloride	13j24-b04-286569	2
Sanacorp	Diltiazem hydrochloride	17337005	18
Sanacorp	Diltiazem hydrochloride	11J28-04511	1
Sanacorp	Diphenhydramine hydro...	12329805	3
Sanacorp	Diphenhydramine hydro...	15072617	2
Sanacorp	Dithranol	APOP-DNL/11/17/001-D	2
Sanacorp	Emulsifying cetostear...	12J29-N04	4
Sanacorp	Eosin Y	12E25-N01	5
Sanacorp	Equisetum stem, cut	13344307	1
Sanacorp	Erythromycin	59/AKD51	1
Sanacorp	Erythromycin	CH.-B.59/AKD50	1
Sanacorp	Erythromycin	16J24-B10	6
Sanacorp	Erythromycin	181101901	1
Sanacorp	Estradiol benzoate	10/AKD24	1
Sanacorp	Estradiol valerate	13264002	1
Sanacorp	Estriol	16317701	5
Sanacorp	European golden rod, ...	12343103	1
Sanacorp	Feverfew herb, cut	12380201	2
Sanacorp	Fructose	K43537721	2
Sanacorp	Gentamicin sulfate	FG1711309	1
Sanacorp	Gentamicin sulfate	13189607	6
Sanacorp	Gentamicin sulfate	17358618	3
Sanacorp	Gentamicin sulfate	18C12-B02-349026	1
Sanacorp	Gentamicin sulfate	161003601	4
Sanacorp	Glucose monohydrate	18101903	11
Sanacorp	Glucose monohydrate	3/AKD100	1
Sanacorp	Glucose monohydrate	16184810	25
Sanacorp	Goat's rue, cut	12017911	1
Sanacorp	Griseofulvin	17L05-B05-345995	1
Sanacorp	Hard fat in pastilles	16274202	2
Sanacorp	Heavy alkaline bismut...	12146909	2
Sanacorp	Hydrochlorothiazide	12A03-B40-26111	1
Sanacorp	Hydrocortisone	1112030-01	9
Sanacorp	Hydrocortisone acetate	18278006	9
Sanacorp	Hydroquinone	12288503	5
Sanacorp	Hydroquinone	16023103	1
Sanacorp	Hyssop herb, cut	11374303	2
Sanacorp	Hyssop herb, cut	13336101	7
Sanacorp	Juniper	12330509	1
Sanacorp	Kelp, cut	11351105	1
Sanacorp	Ladys bedstraw herb,...	11265109	1
Sanacorp	Ladys bedstraw herb,...	12369003	8
Sanacorp	Lidocaine	1603027-02	1
Sanacorp	Lidocaine hydrochlori...	14J07-B04-299434	4
Sanacorp	Lidocaine hydrochlori...	150502302	1
Sanacorp	Lime flower, cut	264045	1
Sanacorp	Liquorice root, peeled	11351906	3
Sanacorp	Lycopodium herb, cut	13101502	2
Sanacorp	Macrogols 4000	15376520	1
Sanacorp	Macrogols 4000	14061707	5
Sanacorp	Mallow leaf, cut	278925	1
Sanacorp	Matricaria flower	12294804	6
Sanacorp	Matricaria flower	12294809	3
Sanacorp	Matricaria flower	1300A120725	2
Sanacorp	Matricaria flower	1290B140123	2
Sanacorp	Meadowsweet herb, cut	4541A120313	2
Sanacorp	Melatonin	APALC/MLT/161101M	1
Sanacorp	Melatonin	APALC/MLT&170301M	1
Sanacorp	Melatonin	APALC/MLT/170301M	2
Sanacorp	Melatonin	MLT/SEP/170901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Melatonin	APMLT170901M	1
Sanacorp	Melatonin	APALC/MLT/170901M	1
Sanacorp	Melatonin	161004002	1
Sanacorp	Melatonin	161004001	2
Sanacorp	Menthol	17328804	9
Sanacorp	Menthol	6/AKD3	1
Sanacorp	Methylcellulose	1717405	1
Sanacorp	Metronidazole	39/AKD18	1
Sanacorp	Metronidazole	16I19-B04-329160	3
Sanacorp	Metronidazole	39/AKD6	1
Sanacorp	Metronidazole	601021826	1
Sanacorp	Metronidazole	39AKD15	1
Sanacorp	Metronidazole	12110422	17
Sanacorp	Metronidazole	CH.-B.39/AKD17	1
Sanacorp	Metronidazole	17C10-B04-338707	1
Sanacorp	Metronidazole	17C10-B04-346280	2
Sanacorp	Miconazole nitrate	CH.-B.4/AKD10	1
Sanacorp	Minoxidil	13179502	7
Sanacorp	Minoxidil	14093702	3
Sanacorp	Minoxidil	120056-01	1
Sanacorp	Minoxidil	14093701	1
Sanacorp	Mometasone furoate, m...	15189801	2
Sanacorp	Mometasone furoate, m...	17D19-B02-344776	1
Sanacorp	Mometasone furoate, m...	MOM001/0317	2
Sanacorp	Mometasone furoate, m...	17D19-B02-339276	1
Sanacorp	Neostigmine bromide	18024902	1
Sanacorp	Nettle leaf, cut	12184707	11
Sanacorp	Nettle leaf, cut	12184003	17
Sanacorp	Nettle root, cut	11265001	2
Sanacorp	Nettle seed	13010602	1
Sanacorp	Omeprazole	17H15-B03-342821	4
Sanacorp	Orris root	11026701	1
Sanacorp	Polymyxin B sulfate	1707013-01	6
Sanacorp	Potassium hydrogen ta...	15354609	2
Sanacorp	Prednicarbate, micron...	17315604	22
Sanacorp	Prednisolone acetate	15398408	11
Sanacorp	Prednisolone acetate	15123-B04-341793	1
Sanacorp	Prednisolone acetate	15123-B04-334688	2
Sanacorp	Prednisolone acetate	18067303	9
Sanacorp	Prednisolone, microni...	1108055-03	13
Sanacorp	Prednisolone, microni...	16H04-B03_332930	1
Sanacorp	Prednisolone, microni...	14/10-B03-302269	2
Sanacorp	Prednisolone, microni...	13j28-b06-289670	1
Sanacorp	Prednisolone, microni...	18AB02	1
Sanacorp	Primrose flower	1541B130808	6
Sanacorp	Progesterone, microni...	15/11-B05-315224	6
Sanacorp	Progesterone, microni...	17D25-B06	1
Sanacorp	Propranolol hydrochlo...	12370712	1
Sanacorp	Propyl-4-hydroxybenzo...	AKD3	1
Sanacorp	Pumpkin seed, peeled	16081212	1
Sanacorp	Raspberry leaf, cut	13038001	6
Sanacorp	Raspberry leaf, cut	12362703	9
Sanacorp	Resocinol	16D29-Bo1-323046	1
Sanacorp	Ribwort plantain, cut	263401	1
Sanacorp	Roselle	12267508	2
Sanacorp	Salicylic acid	12134707	7
Sanacorp	Salicylic acid	17C21-B09-339446	1
Sanacorp	Salicylic acid	68/AKD2	1
Sanacorp	Salicylic acid	123036223	1
Sanacorp	Senna leaf, cut	13136701	4

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Sildenafil citrate	16049811	3
Sanacorp	Sildenafil citrate	14H11-B01-297519	2
Sanacorp	Sodium acetate trihyd...	13295601	2
Sanacorp	Sodium citrate	1436407	1
Sanacorp	Sodium dodecyl sulfate	5/AKD1	1
Sanacorp	Sodium dodecyl sulfate	5/AKD6	1
Sanacorp	Sodium edetate	1505603	1
Sanacorp	Sodium edetate	12G25-N12	5
Sanacorp	Sodium edetate	6/AKD14	1
Sanacorp	Soy lecithin (instant...	14196305	1
Sanacorp	Soy lecithin (instant...	17322405	1
Sanacorp	Spironolactone	18208004	1
Sanacorp	Spironolactone	12049130	5
Sanacorp	Spironolactone	18004913	4
Sanacorp	Spironolactone	1212034-01	1
Sanacorp	Spironolactone	11042612	1
Sanacorp	Spironolactone	12049103	3
Sanacorp	Sucrose	AKD5	1
Sanacorp	Sucrose	14340005	2
Sanacorp	Tannin	CH.-B.AKD60	1
Sanacorp	Tannin	12336011	2
Sanacorp	Tannin	12171102	1
Sanacorp	Tetracaine hydrochlor...	17G03-B07-341853	1
Sanacorp	Tetracaine hydrochlor...	15/AKD1	1
Sanacorp	Tetracycline hydrochl...	100403401	1
Sanacorp	Thymol	9E14-N16	1
Sanacorp	Thymol	1135816	1
Sanacorp	Tormentil	12040712	3
Sanacorp	Tosylchloramide sodium	12H09-N01	2
Sanacorp	Tretinoin	1411016-01	11
Sanacorp	Tretinoin	1411016-02	12
Sanacorp	Tretinoin	12/06-N01	1
Sanacorp	Tretinoin	12320302	1
Sanacorp	Triamcinolone acetone...	NY0485M	2
Sanacorp	Triamcinolone acetone...	17H03-B12-345553	2
Sanacorp	Triamcinolone acetone...	17358305	16
Sanacorp	Triamcinolone acetone...	13124807	4
Sanacorp	Triamcinolone acetone...	1326N11	1
Sanacorp	Triamcinolone acetone...	15430108	9
Sanacorp	Triamcinolone, micron...	15106907	1
Sanacorp	Triamcinolone, micron...	15106916	1
Sanacorp	Triclosan	28/AKD4	1
Sanacorp	Triclosan	13391807	6
Sanacorp	Triclosan	CH.-B.22/AKD25	1
Sanacorp	Triclosan	26/AKD1	1
Sanacorp	Trometamol	3686G018	1
Sanacorp	Trometamol	8386H037	1
Sanacorp	Trometamol	17320804	12
Sanacorp	Trometamol	12113506	9
Sanacorp	Urea	62/AKD50s	1
Sanacorp	Urea	1212007	1
Sanacorp	Urea	63/AKD1	1
Sanacorp	Urea	13064212	27
Sanacorp	Urea	1405028-01	1
Sanacorp	Valerian root, cut	12033804	9
Sanacorp	Vanillin	12H09-N06	4
Sanacorp	Wax, yellow	12/AKD2	1
Sanacorp	Wax, yellow	13378619	14
Sanacorp	Wax, yellow	12257009	8
Sanacorp	Wild pansy, cut	12221-01694	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Sanacorp	Wormwood herb, cut	12103002	2
Sanacorp	Xylometazoline hydroc...	15037906	2
SANACORP	Ascorbic acid	201401318	1
Sanacorp 28.06.2016	Sodium edetate	15D01-B05-315877	1
Sanacorp 27.10.15	Metronidazole	15A20-B07-303298	10
Sanacorp WE: 12.08.17...	Calcium citrate tetra...	711105	1
Sanacorp WE: 28.09.17...	Gentamicin sulfate	15258825	2
Sanacorp, 01.03.17, 2...	Prednicarbate, micron...	1687410	1
Sanacorp, 39,57EUR, 27...	Alfatradiol	700670-0011	1
Sanacorp/ euRho	Potassium citrate	1204016-01	2
Sanacorp/ F	Griseofulvin	12B06-N07	1
Sanacorp/Phoenix	Spironolactone	13124315	7
Sanacorp;19.03.15;87,...	Clobetasol propionate	14C03-B07-293563	2
sanacorp;3,32;18.06.15	Coriander	13299495	1
Sanacorp04	Clioquinol	1612020-02	3
Sanacorpaelo /Sanacorp	Cetyl alcohol	1AKS11	1
Sanofi Chimie/ Gehe	Dexamethasone acetate	18173004	3
schubert / Phoenix	Lactose monohydrate	1334PLV-1775	1
schubert / Phoenix	Tartaric acid	1014PLV-755	1
Shandong northwest ma...	Alfatradiol	AEC151211	1
Shandong Northwest Ma...	Alfatradiol	AEC160522	1
Shandong Northwest Ma...	Estriol	ETL160419	1
Shandong Northwest Ma...	Estriol	ETL160302	1
Shandong Northwest Ma...	Finasteride (form III)	FNA160604	1
Shandong Northwest Ma...	Testosterone	GTJ160510	1
Shop	Calcium lactate penta...	1111001829	1
Sidonia	Bismuth gallath, basic	17004107	6
Sidonia	Bismuth gallath, basic	15J-26-B02-325450	1
Sidonia	Metoprolol tartrate	405005RIIC/09062017	1
Sidonia	Prednisolon dihydroge...	17I27-B03-342560	1
Sidonia	Prednisolon dihydroge...	17I27-B03-344246	2
Sidonia	Prednisolon dihydroge...	17I27-B03-345405	7
siehe Protokoll	Aniseed	sieheProtokoll	1
siehe Protokoll	Cinnamon, cut	sieheProtokoll	1
siehe Protokoll	Diclofenac sodium	15169617	1
siehe Protokoll (Galke)	Bitter fennel	sieheProtokoll	2
siehe WE-Prüfprotokoll	8-quinolinol sulfate ...	96/09	1
siehe WE-Prüfprotokoll	Bismuth gallath, basic	4/09	1
siehe WE-Prüfprotokoll	Neomycin sulfate	19/12	1
siehe WE-Prüfprotokoll	Potassium sorbate	61/11	1
siehe WE-Prüfprotokoll	Procaine hydrochloride	39/12	1
siehe WE-Prüfprotokoll	Propyl-4-hydroxybenzo...	84/11	1
siehe WE-Prüfprotokoll	Soy lecithin (instant...	31/14	1
siehe WE-Prüfprotokoll	Uridine	54/11	1
Sigma Adrich/Alliance...	Metoprolol tartrate	PHR1076	1
Sigma-Aldrich	Sodium butyrate	Lot.091M-5302V	1
Sigma-Aldrich	Sodium proprionate	SLBK5249V	1
Sigma-Aldrich GmbH	Metoprolol tartrate	LRAA6905	4
Sino Phyto	Agrimony, cut	12382905	5
Soldan	Sucrose	25.03.2014	1
Soldan	Sucrose	232027	1
Soldan	Sucrose	243983	1
Sonnen-Apotheke Wefen...	Glucose	17/16	1
Sonnen-Apotheke Wefen...	Glucose	34/16	1
Sonnenberg-Apotheke	Captopril	14B26-B01-291557	3
Sonnenberg-Apotheke	Glucose monohydrate	140500401	3
Sonnenberg-Apotheke	Potassium sorbate	13224604	18
Sonnenberg-Apotheke	Prednicarbate, micron...	1607012-01	9
Spagyra	Sucrose	5030921072	1
Spagyra	Sucrose	1905	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Spagyra	Sucrose	1691	1
Spagyra	Sucrose	17072120536	1
Spagyra	Sucrose	1918	1
Spagyra	Sucrose	17072120537	1
Spagyra	Sucrose	17070420968	1
Spagyra	Sucrose	17082320366	1
Spagyra	Sucrose	5030921078	1
Spagyra	Sucrose	2007	1
Spagyra	Sucrose	5030921056	1
Spagyra	Sucrose	5030920922	1
Spagyra	Sucrose	5030921081	1
Spagyra	Sucrose	5030921067	1
Spagyra	Sucrose	16110320884	1
Spagyra	Sucrose	2923	1
Spagyra	Sucrose	5030920895	1
Spagyra	Sucrose	17072120538	1
Spagyra	Sucrose	5030921080	1
Spagyra	Sucrose	5030921152	1
Spagyra	Sucrose	17090620617	1
Spagyra	Sucrose	5030921094	1
Spagyra	Sucrose	5030921090	1
Spagyra	Sucrose	5030920908	1
Spagyra	Sucrose	10032921256	1
Spagyra	Sucrose	1911	1
Spagyra	Sucrose	5030921060	1
Spagyra	Sucrose	16083120745	1
Spagyra	Sucrose	17080920574	1
Spagyra	Sucrose	5030920924	1
Spagyra	Sucrose	5030921065	1
Spagyra	Sucrose	5030920893	1
Spagyra	Sucrose	15032620414	1
Spagyra	Sucrose	16102120418	1
Spagyra	Sucrose	5030921122	1
Spagyra	Sucrose	5030921020	1
Spagyra	Sucrose	13061020526	1
Spagyra	Sucrose	5030921019	1
Spagyra	Sucrose	13061020529	1
Spagyra	Sucrose	13061020524	1
Spagyra	Sucrose	5030921024	1
Spagyra	Sucrose	13061020528	1
Spagyra	Sucrose	13061220385	2
Spagyra	Sucrose	13050621773	2
Spagyra	Sucrose	13070420937	2
Spagyra	Sucrose	5030921038	1
Spagyra	Sucrose	5030921082	1
Spagyra	Sucrose	5030921010	1
Spagyra	Sucrose	13070420936	2
Spagyra	Sucrose	13040921431	1
Spagyra	Sucrose	5030921039	1
Spagyra	Sucrose	13072620834	1
Spagyra	Sucrose	16032920732	1
Spagyra	Sucrose	16042522126	1
Spagyra	Sucrose	10073020385	1
Spagyra	Sucrose	13061020525	1
Spagyra	Sucrose	5030921042	1
Spagyra	Sucrose	5030920999	1
Spagyra	Sucrose	5030921026	1
Spagyra	Sucrose	12103020975	1
Spagyra	Sucrose	15082720411	1
Spagyra	Sucrose	13040921430	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Spagyra	Sucrose	5030921031	1
Spagyra	Sucrose	5030921030	1
Spagyra	Sucrose	13040921424	1
Spagyra	Sucrose	13012920302	1
Spagyra	Sucrose	13020620240	1
Spagyra	Sucrose	13030621373	1
Spagyra	Sucrose	5030921014	1
Spagyra	Sucrose	5030921006	1
Spagyra	Sucrose	13041920154	1
Spagyra	Sucrose	13041920152	2
Spagyra	Sucrose	5030921032	1
Spagyra	Sucrose	5030921041	1
Spagyra	Sucrose	5030921036	1
Spagyra	Sucrose	13051520939	1
Spagyra	Sucrose	5030921040	1
Spagyra	Sucrose	13051520941	1
Spagyra	Sucrose	5030921028	1
Spagyra	Sucrose	13040921427	1
Spagyra	Sucrose	5030921064	1
Spagyra	Sucrose	13070420938	1
Spagyra	Sucrose	16053020827	1
Spagyra	Sucrose	5030920989	1
Spagyra	Sucrose	5030920995	1
Spagyra	Sucrose	5030920881	1
Spagyra	Sucrose	5030921043	1
Spagyra	Sucrose	5030920883	1
Spagyra	Sucrose	5030921053	1
Spagyra	Sucrose	5030921101	1
Spagyra	Sucrose	5030921002	1
Spagyra	Sucrose	5030921100	1
Spagyra	Sucrose	5030921221	1
Spagyra	Sucrose	5030921096	1
Spagyra	Sucrose	5030921012	1
Spagyra	Sucrose	5030921015	1
Spagyra	Sucrose	5030921016	1
Spagyra	Sucrose	5030920923	1
Spagyra	Sucrose	3750	1
Spagyra	Sucrose	18081720533	1
Spagyra	Sucrose	5030921054	1
Spagyra	Sucrose	5030921004	1
Spagyra	Sucrose	5030921037	1
Spagyra	Sucrose	16051820801	1
Spagyra	Sucrose	16050220821	1
Spagyra	Sucrose	16050421023	1
Spagyra	Sucrose	5030921075	1
Spagyra	Sucrose	5030920873	1
Spagyra	Sucrose	7062220875	1
Spagyra	Sucrose	5030920997	1
Spagyra	Sucrose	5030920998	1
Spagyra	Sucrose	5030920885	1
Spagyra	Sucrose	16070421669	1
Spagyra	Sucrose	12082320717	1
Spagyra	Sucrose	5030920996	1
Spagyra Österreich	Sucrose	15112321097	1
Spagyra/Österreich	Sucrose	17021620307	1
Spagyra/Österreich	Sucrose	H18-0021	1
Spagyra/Österreich	Sucrose	H18-0210	1
Spagyra/Österreich	Sucrose	17020820816	1
Spagyra/Österreich	Sucrose	1809272503	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Spagyra/Österreich	Sucrose	3183	1
Spagyra/Österreich	Sucrose	18021921250	1
Spagyra/Österreich	Sucrose	1721420709	1
Spagyra/Österreich	Sucrose	17032120927	1
Spagyra/Österreich	Sucrose	18092720503	1
Spagyra/Österreich	Sucrose	17070320627	1
Spagyra/Österreich	Sucrose	17021620309	1
Spagyra/Österreich	Sucrose	17041821026	1
Spagyra/Österreich	Sucrose	17022022188	1
Spagyra/Österreich	Sucrose	3303	1
Spagyra/Österreich	Sucrose	2382	1
Spagyra/Österreich	Sucrose	17092620773	1
Spagyra/Österreich	Sucrose	17020820810	1
Spagyra/Österreich	Sucrose	2307	2
Spagyra/Österreich	Sucrose	17041821027	1
Spagyra/Österreich	Sucrose	17121420710	1
Spagyra/Österreich	Sucrose	17091520139	1
Spagyra/Österreich	Sucrose	17032820797	1
Spagyra/Österreich	Sucrose	17032321132	1
Spagyra/Österreich	Sucrose	16122820921	1
Spagyra/Österreich	Sucrose	2516	1
Spagyra/Österreich	Sucrose	2963	1
Spagyra/Österreich	Sucrose	H17-0047	1
Spagyra/Österreich	Sucrose	3082	1
Spagyra/Österreich	Sucrose	14012420377	1
Spagyra/Österreich	Sucrose	18072420037	1
Spagyra/Österreich	Sucrose	1624	1
Spagyra/Österreich	Sucrose	18041120712	1
Spagyra/Österreich	Sucrose	14050621091	1
Spagyra/Österreich	Sucrose	13091920456	1
Spagyra/Österreich	Sucrose	18090520789	1
Spagyra/Österreich	Sucrose	18051421071	1
Spagyra/Österreich	Sucrose	18051421069	1
Spagyra/Österreich	Sucrose	12112721050	1
Spagyra/Österreich	Sucrose	14010721509	1
Spagyra/Österreich	Sucrose	14060353342	1
Spagyra/Österreich	Sucrose	13050621774	1
Spagyra/Österreich	Sucrose	18042720542	1
Spagyra/Österreich	Sucrose	14071421137	2
Spagyra/Österreich	Sucrose	18042720555	2
Spagyra/Österreich	Sucrose	13021920149	1
Spagyra/Österreich	Sucrose	3071	1
Spagyra/Österreich	Sucrose	13102520353	1
Spagyra/Österreich	Sucrose	18041180714	1
Spagyra/Österreich	Sucrose	18091020640	1
Spagyra/Österreich	Sucrose	14041421810	1
Spagyra/Österreich	Sucrose	14060353343	1
Spagyra/Österreich	Sucrose	18032320386	1
Spagyra/Österreich	Sucrose	18090520787	3
Spagyra/Österreich	Sucrose	18090520788	1
Spagyra/Österreich	Sucrose	18032320388	1
Spagyra/Österreich	Sucrose	5030920991	1
Spagyra/Österreich	Sucrose	18032320385	1
Spagyra/Österreich	Sucrose	5030921018	1
Spagyra/Österreich	Sucrose	5030921008	1
Spagyra/Österreich	Sucrose	5030920993	1
Spagyra/Österreich	Sucrose	14070721041	1
Spagyra/Österreich	Sucrose	10070521044	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Spagyra/Österreich	Sucrose	18032320387	1
Spangro/gepepharm	Bitter fennel	15007602	2
Spangropharm	Birch leaf, cut	264603	1
Spangropharm	Glycerol monostearate...	80402468	1
Spangropharm	Lactose monohydrate	91444379	1
Spangropharm	Matricaria flower	11169606	1
Spangropharm	Metronidazole	11319818	4
Spangropharm	Mistletoe, cut	10273408	1
Spangropharm	Tartaric acid	10073111	1
spinnrad	Potassium sorbate	1118331022	1
spinnrad	Potassium sorbate	2233102	1
Spinnrad	Black cumin seed	11931059	1
Spinnrad	Cocoa butter	512300202	1
Spinnrad	Cocoa butter	1017300203	1
Spinnrad	Potassium sorbate	1117331022	5
Spinnrad	Potassium sorbate	312331021	2
Spinnrad	Potassium sorbate	317331021	5
Spinnrad	Potassium sorbate	815331022	2
Spinnrad	Potassium sorbate	214331021	2
Spinnrad	Wax, bleached	317048061	1
Spinnrad	Wax, bleached	118048061	1
Spinnrad GmbH	Potassium sorbate	111033102	3
Spinnrad GmbH	Potassium sorbate	912331023	7
Spinnrad GmbH	Potassium sorbate	1114331022	6
Spinnrad/Gehe	Potassium sorbate	315331021	4
Spinnrad/Phoenix	Black cumin seed	613310592	1
Spinnrad/Sannacorp	Potassium sorbate	616331022	4
Spinrad/Phönix	Potassium sorbate	218331021	1
Spirig	Testosterone	13406003	1
St. Erhard Apotheke	Sucrose	18090520790	1
St. Erhard Apotheke	Sucrose	18090520795	1
St. Erhard Apotheke	Sucrose	1809052072	1
St. Erhard Apotheke	Sucrose	18090520792	1
St. Erhard Apotheke	Sucrose	18090520794	1
St. Erhard Apotheke	Sucrose	18090520791	1
St. Erhard Apotheke	Sucrose	18090520793	1
Staufenpharma	Sucrose	3314E-01586	1
Staufen-Pharma	Sucrose	161112	2
stumpf	Magnesium chloride he...	1301c012-277804	1
stumpf	Mannitol	1104028-02	1
Stumpf	Povidone iodine	11168612	2
Supriya LifeScience L...	Pheniramine maleate	SLL/P/1213063	1
Supriya LifeScience L...	Pheniramine maleate	14C25-B01	1
Swati/Noweda	Lidocaine hydrochlori...	LDH/116007	1
Symbiotica Spec.	Betamethasone valerate	1208014-04	9
Symbiotica Spec./Phönix	Clobetasol propionate	CL0-B12153/0711	1
Symbiotica Speciality	Testosterone propionate	15444702	8
Synop WE:18.12.17 E...	L-isoleucine	17H14-B04-339009	1
Synopharm	Fructose	202A052	1
Synopharm/Anzag	Tannin	3411I-01052	1
Synopharm/Anzag	Tetracaine hydrochlor...	2712Q-02039	1
Taoasis	Ascorbic acid	810049	1
Taoasis	Tosylchloramide sodium	1623607	1
Taoasis/Phoenix	Neomycin sulfate	1502004-01	6
test	Urea	test	1
Test	Sodium citrate		7
Th. Geyer	Dequalinium chloride	15B16-B04-354618	1
Th. Geyer	Lipoic acid (a-) (H10...	15001526/0	1
Th. Geyer	Metoprolol tartrate	56392-17-7	1
THC Pharm	Cannabidiol	CBDAP1901	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
THC Pharm	Cannabidiol	H150.2318-1	1
THC Pharm	Cannabidiol	H150218-1	1
THC Pharm	Cannabidiol	CBD-API16-002	1
THC Pharm	Cannabidiol	H150.218-3	5
THC Pharm	Cannabidiol	CBD-API16-004	6
THC Pharm GmbH	Cannabidiol	CBD-API1802	5
THC Pharm GmbH	Cannabidiol	H150.218-1	9
THC Pharm GmbH	Cannabidiol	H160.718-1	35
THC Pharma	Cannabidiol	CBDAP/1801	1
THC Pharma	Cannabidiol	CBD-API17-005	12
THC Pharma	Cannabidiol	CBD-API17-002	10
THC Pharma	Cannabidiol	H150.218-4	4
THC Pharma	Cannabidiol	H150.215-4	1
THC Pharma	Cannabidiol	H150.218.-2	1
THC Pharma	Cannabidiol	H150.218-2	5
THC Pharma	Cannabidiol	CBD-API17-004	15
THC Pharma	Cannabidiol	H.180.119-2	1
THC Pharma (direkt)	Cannabidiol	CBD-API1707	4
THC Pharma (direkt)	Cannabidiol	CBD-API1708	2
THC/ Direkt	Cannabidiol	CBD-API1706	5
THC-Pharm	Cannabidiol	cbd-adi17-002	1
THC-Pharm	Cannabidiol	110.216-2	8
THC-Pharm	Cannabidiol	V040.316-3	3
THC-Pharm	Cannabidiol	V110.417-1	1
THC-Pharm	Cannabidiol	CBDAP/17-003	3
THC-Pharma	Cannabidiol	CBD-API17-003	6
THC-Pharma	Cannabidiol	H041.018-3	28
THC-Pharma	Cannabidiol	H150,218-2	1
THC-Pharma	Cannabidiol	CBDAP/1708	1
un	Vitamin E NAT (H1107)	25022	5
unbekannt	Magnesium chloride he...	unbekannt	1
unisan	Aminomix NAC (H1002)	15000297	1
unisan	Aminomix NAC (H1002)		1
unisan	Antioxidants NAT comp...	15000877	2
unisan	Chromium (H1012)		2
unisan	Inositol hexanicotina...	15000013	2
unisan	Lipoic acid (a-) (H10...	66601	4
unisan	Lipoic acid (a-) (H10...	H100031	1
unisan	Magnesium (H1060)	43702	2
unisan	Taurine (H1084)	14001767	2
unisan	Tryptophan (L-) (H1053)		1
unisan	Vitamin C (H1104)	14001794	1
Unisan	Acacia gum (H1001)	16000363/0	1
Unisan	Aminomix NAC (H1002)	17000766	3
Unisan	Aminomix NAC (H1002)	165001017/0	1
Unisan	Aminomix NAC (H1002)	971701	1
Unisan	Aminomix NAC (H1002)	19101	1
Unisan	Aminomix NAC (H1002)	14000258	3
Unisan	Aminomix NAC (H1002)	15000298/0	1
Unisan	Aminomix NAC (H1002)	140000259	1
Unisan	Aminomix NAC (H1002)	14000259	1
Unisan	Aminomix NAC (H1002)	15000320/0	3
Unisan	Aminomix NAC (H1002)	15000295/0	3
Unisan	Aminomix NAC (H1002)	41401	3
Unisan	Antioxidants NAT comp...	14000961	2
Unisan	Antioxidants NAT comp...	67501	2
Unisan	Antioxidants NAT comp...	1400042	1
Unisan	Antioxidants NAT comp...	15000877/0	1
Unisan	Antioxidants NAT comp...	11001	1
Unisan	Antioxidants NAT comp...	14000044	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	Antioxidants NAT comp...	30903	3
Unisan	Antioxidants NAT comp...	1034030904	2
Unisan	Antioxidants NAT comp...	1851	1
Unisan	Antioxidants NAT comp...	17000264	3
Unisan	Antioxidants NAT comp...	17000260	3
Unisan	Antioxidants NAT comp...	30904	1
Unisan	Antioxidants NAT comp...	16000970/0	4
Unisan	Antioxidants NAT comp...	19901	1
Unisan	Antioxidants NAT comp...	51901	3
Unisan	Antioxidants NAT comp...	14000963	5
Unisan	Antioxidants NAT comp...	16000968/0	1
Unisan	Antioxidants NAT comp...	67201	5
Unisan	Antioxidants NAT comp...	67601	2
Unisan	Antioxidants NAT comp...	15000880	1
Unisan	Antioxidants NAT comp...	14001851/0	1
Unisan	Antioxidants NAT comp...	14001848/0	1
Unisan	Antioxidants NAT comp...	51902	2
Unisan	Antioxidants NAT comp...	14001847/1	1
Unisan	Antioxidants NAT comp...	16000378	1
Unisan	Antioxidants NAT comp...	16000966/0	2
Unisan	Antioxidants NAT comp...	16000966	1
Unisan	Antioxidants NAT comp...	14001847/0	2
Unisan	Antioxidants NAT comp...	68201	4
Unisan	Antioxidants NAT comp...	52001	2
Unisan	Antioxidants NAT comp...		1
Unisan	Antioxidants NAT comp...	19901	2
Unisan	Arginine (L-) (H1036)	66501	7
Unisan	Arginine (L-) (H1036)	27802	1
Unisan	Arginine (L-) (H1036)	44301	2
Unisan	Arginine (L-) (H1036)	14000563	1
Unisan	Arginine (L-) (H1036)	14000564/0	1
Unisan	Arginine (L-) (H1036)	14000564/2	1
Unisan	Arginine (L-) (H1036)	16000949/0	1
Unisan	Arginine (L-) (H1036)	17000877	3
Unisan	Arginine (L-) (H1036)	1324950601	1
Unisan	Betamethasone, micron...	Retest09K25-N13	1
Unisan	Bioflavonoid complex ...	1281683901	1
Unisan	Bioflavonoid complex ...	14001675/4	1
Unisan	Bioflavonoid complex ...	14001675	2
Unisan	Bioflavonoid complex ...	Pr.Nr.1452	1
Unisan	Bioflavonoid complex ...	1284648802	1
Unisan	Bioflavonoid complex ...	42501	1
Unisan	Biotin (H1007)	14000862	2
Unisan	Biotin (H1007)	15000735/5	1
Unisan	Biotin (H1007)	983601	1
Unisan	Biotin (H1007)	16001396	3
Unisan	Biotin (H1007)	29603	6
Unisan	Biotin (H1007)	16001396/0	1
Unisan	Biotin (H1007)	15000735/0	1
Unisan	Biotin (H1007)	15000735	2
Unisan	Biotin (H1007)	14000862/2	2
Unisan	Calcium (H1008)	8901	8
Unisan	Calcium (H1008)	55101	3
Unisan	Calcium (H1008)	17000201	1
Unisan	Calcium (H1008)	956001	1
Unisan	Calcium (H1008)	15001666/0	1
Unisan	Calcium (H1008)	15001666	1
Unisan	Calcium (H1008)	14001938/0	1
Unisan	Calcium (H1008)	14001934/0	1
Unisan	Calcium (H1008)	14000108	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	Calcium (H1008)	14000110	2
Unisan	Calcium (H1008)	55301	3
Unisan	Calcium (H1008)	15001667	2
Unisan	Calcium (H1008)	1105008901	1
Unisan	Calcium (H1008)	15001672/0	1
Unisan	Calcium (H1008)	55601	2
Unisan	Carnitine (L-) (H1038)	15000053/3	1
Unisan	Carnitine (L-) (H1038)	15000054/0	1
Unisan	Carnitine (L-) (H1038)	15000054	1
Unisan	Carnitine (L-) (H1038)	150000540	1
Unisan	Carnitine (L-) (H1038)	15000053	1
Unisan	Carnitine (L-) (H1038)	15000053/0	2
Unisan	Carnitine (L-) (H1038)	16000450	2
Unisan	Carnitine (L-) (H1038)	15000052	1
Unisan	Carnitine (L-) (H1038)	26010	1
Unisan	Carnitine (L-) (H1038)	130605301	2
Unisan	Carnitine (L-) (H1038)	130605300	2
Unisan	Carnitine (L-) (H1038)	1184026020	3
Unisan	Carnitine (L-) (H1038)	130326002	1
Unisan	Carnitine (L-) (H1038)	965901	1
Unisan	Carnitine (L-) (H1038)	130605300/3	1
Unisan	Carnitine (L-) (H1038)	16000451	2
Unisan	Carnitine (L-) (H1038)		2
Unisan	Catechin extract (H10...	17000409	2
Unisan	Catechin extract (H10...	1571025702	2
Unisan	Catechin extract (H10...	1571681501	1
Unisan	Catechin extract (H10...	25702	1
Unisan	Catechin extract (H10...	130208020	2
Unisan	Catechin extract (H10...	681501	1
Unisan	Catechin extract (H10...	62801	1
Unisan	Chaga (H1011)	14001688	1
Unisan	Chaga (H1011)	64001	2
Unisan	Chaga (H1011)	16000235	1
Unisan	Chondroitin sulfate ...	1820001801	2
Unisan	Chondroitin sulfate ...	1801	1
Unisan	Chondroitin sulfate ...	1801	1
Unisan	Chromium (H1012)	16000520/0	1
Unisan	Chromium (H1012)	977201	10
Unisan	Chromium (H1012)	1582977201	4
Unisan	Chromium (H1012)	1581977201	1
Unisan	Clioquinol	14/11-B03	1
Unisan	Coenzyme Q10 (H1013)	992501	1
Unisan	Coenzyme Q10 (H1013)	16000402	1
Unisan	Coenzyme Q10 (H1013)	16000402/0	3
Unisan	Coenzyme Q10 (H1013)	H101311	1
Unisan	Coenzyme Q10 (H1013)	16001598	1
Unisan	Coenzyme Q10 (H1013)	15601	9
Unisan	Coenzyme Q10 (H1013)	14001546/1	2
Unisan	Copper (H1032)	15000211	1
Unisan	Copper (H1032)	7901	1
Unisan	Copper (H1032)	46401	3
Unisan	Copper (H1032)	15000211/0	1
Unisan	Creatine (H1031)	15000750	1
Unisan	Creatine (H1031)	25003	1
Unisan	Creatine (H1031)	25003	1
Unisan	Creatine (H1031)	15000750/0	1
Unisan	Curcumin and piperine...	14001017/6	1
Unisan	Curcumin and piperine...	14001017	1
Unisan	Curcumin and piperine...	72001	6
Unisan	Curcumin and piperine...	17000259	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	D-galactose	432619PF5	1
Unisan	D-galactose	382619PF5	1
Unisan	D-galactose	442619PF5	1
Unisan	D-ribose	JXCZSJ13D-013	2
Unisan	Folic acid (H1016)	30201	9
Unisan	Folic acid (H1016)	1613977701	2
Unisan	Folic acid (H1016)	H101621	1
Unisan	Folic acid (H1016)	14000639/1	2
Unisan	Folic acid (H1016)	977701	2
Unisan	Folic acid (H1016)	16000931	3
Unisan	Folic acid (H1016)	14000638	3
Unisan	Folic acid (H1016)	18000337	2
Unisan	Folic acid (H1016)	16001098	4
Unisan	Garcinia cambogia (H1...	17000873	1
Unisan	Garcinia cambogia (H1...	56601	1
Unisan	Garcinia cambogia (H1...	1221007401	2
Unisan	Garcinia cambogia (H1...	16001252/0	1
Unisan	Ginseng (H1018)	14002047/2	1
Unisan	Ginseng (H1018)	17000285	2
Unisan	Ginseng (H1018)	14002047/1	1
Unisan	Ginseng (H1018)	67801	3
Unisan	Glucosamine sulfate (...)	75401	1
Unisan	Glucosamine sulfate (...)	16000948	1
Unisan	Glucosamine sulfate (...)	17000073/0	1
Unisan	Glucosamine sulfate (...)	1823954901	2
Unisan	Glutamine (L-) (H1043)	15000018/0	1
Unisan	Glutamine (L-) (H1043)	15000018	1
Unisan	Glutamine (L-) (H1043)	16000940/0	2
Unisan	Glutamine (L-) (H1043)	16000940	1
Unisan	Glutamine (L-) (H1043)	16000941	3
Unisan	Glutamine (L-) (H1043)	16000942	4
Unisan	Glutamine (L-) (H1043)		1
Unisan	Glutamine (L-) (H1043)	15000017	1
Unisan	Glutamine (L-) (H1043)	15000017/0	2
Unisan	Glutamine (L-) (H1043)	1814026901	5
Unisan	Glutamine (L-) (H1043)	61501/4	1
Unisan	Glutamine (L-) (H1043)	26903	2
Unisan	Glutamine (L-) (H1043)	61501	3
Unisan	Glutamine (L-) (H1043)	61501	7
Unisan	Glutathione (H1020)	16000938	1
Unisan	Glutathione (H1020)	25703	1
Unisan	Glutathione (H1020)	40201	4
Unisan	Glycine (H1021)	14001730	1
Unisan	Glycine (H1021)	16000933	2
Unisan	Inositol hexanicotina...	14000920	1
Unisan	Iron (H1015)	14000285/5	1
Unisan	Iron (H1015)	1114603206	1
Unisan	Iron (H1015)	40901	11
Unisan	Iron (H1015)	14000285/1	2
Unisan	Iron (H1015)	15001259	2
Unisan	Iron (H1015)	15001258/2	2
Unisan	Iron (H1015)	15001258/0	3
Unisan	Iron (H1015)	998701	6
Unisan	Iron (H1015)	40901	2
Unisan	Isoflavones (H1025)	15000205/3	1
Unisan	Isoflavones (H1025)	15000205	1
Unisan	Isoflavones (H1025)	15000205/0	3
Unisan	Isoleucine (L-) (H1045)	1805968701	1
Unisan	Isoleucine (L-) (H1045)	27701	1
Unisan	Isoleucine (L-) (H1045)	66301	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	Isoleucine (L-) (H1045)	15001640	1
Unisan	L-5-HTP (H1035)	26102	1
Unisan	L-5-HTP (H1035)	25006	1
Unisan	L-5-HTP (H1035)	25006	2
Unisan	L-5-HTP (H1035)	25006/12	1
Unisan	L-5-HTP (H1035)	15001503/5	2
Unisan	Leucine (L-) (H1046)	15000283/0	1
Unisan	Leucine (L-) (H1046)	15000283	2
Unisan	Lipoic acid (a-) (H10...	993101	2
Unisan	Lipoic acid (a-) (H10...	14000497	1
Unisan	Lipoic acid (a-) (H10...	25701	10
Unisan	Lipoic acid (a-) (H10...	1312993101	1
Unisan	Lycopene (H1058)	4501	4
Unisan	Lycopene (H1058)	53501/7	1
Unisan	Lycopene (H1058)	53501	8
Unisan	Lycopene (H1058)	H105811	1
Unisan	Lycopene (H1058)	14000803	1
Unisan	Lycopene (H1058)	47301	4
Unisan	Lysine (L-) (H1047)	15000372/0	1
Unisan	Lysine (L-) (H1047)	14000944	1
Unisan	Lysine (L-) (H1047)	1800017302	1
Unisan	Lysine (L-) (H1047)	65401/4	1
Unisan	Magnesium (H1060)	H106021	1
Unisan	Magnesium (H1060)	69801	1
Unisan	Magnesium (H1060)	69701	6
Unisan	Magnesium (H1060)	16000272	3
Unisan	Magnesium (H1060)	69801	3
Unisan	Magnesium (H1060)	14001895	6
Unisan	Magnesium (H1060)	69801/2	1
Unisan	Magnesium (H1060)	14001893	1
Unisan	Magnesium (H1060)	973501	3
Unisan	Magnesium (H1060)	43601	4
Unisan	Magnesium (H1060)	16000265/0	1
Unisan	Magnesium (H1060)	16000265	1
Unisan	Magnesium (H1060)	16000265/1	1
Unisan	Magnesium (H1060)	43501	3
Unisan	Methylsulfonylmethane...	25013	1
Unisan	Methylsulfonylmethane...	979401	1
Unisan	Methylsulfonylmethane...	25013	1
Unisan	Methylsulfonylmethane...	15001371	2
Unisan	Methylsulfonylmethane...	16000449/1	1
Unisan	Methylsulfonylmethane...	26101	1
Unisan	Methylsulfonylmethane...	56201	2
Unisan	OPC Grape seed (H1066)	16001439/0	1
Unisan	OPC Grape seed (H1066)		1
Unisan	OPC Grape seed (H1066)	68301	4
Unisan	OPC Grape seed (H1066)	18401	2
Unisan	Ornithine (L-) (H1048)	16000935/0	3
Unisan	Ornithine (L-) (H1048)	979301	2
Unisan	Ornithine (L-) (H1048)	1816979301	1
Unisan	Paraffin, hard	11166911	3
Unisan	Phenylalanine (L-) (H...	74401/0	1
Unisan	Phenylalanine (L-) (H...	25009	7
Unisan	Phenylalanine (L-) (H...	16001096	2
Unisan	Potassium (H1030)	15000597/0	2
Unisan	Potassium (H1030)	39301	10
Unisan	Potassium (H1030)	944201	3
Unisan	Potassium (H1030)	15000600	2
Unisan	Potassium (H1030)	15000598	2
Unisan	Potassium (H1030)	15000599	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	Potassium (H1030)	288903	1
Unisan	Potassium (H1030)	15000597/1	1
Unisan	Potassium (H1030)	15000599/0	1
Unisan	Potassium (H1030)	69401	1
Unisan	Potassium (H1030)	69301	3
Unisan	Potassium (H1030)	130328001	3
Unisan	Potassium (H1030)	69201	3
Unisan	Potassium (H1030)	28903	1
Unisan	Potassium (H1030)	69401/0	1
Unisan	Proline (-L) (H1051)	66401	2
Unisan	Rhodiola (H1064)	15001045/8	1
Unisan	Selenium (H1067)	16001349	1
Unisan	Selenium (H1067)	1144010201	5
Unisan	Selenium (H1067)	10201	20
Unisan	Selenium (H1067)	1144933301	1
Unisan	Selenium (H1067)	65701	5
Unisan	Selenium (H1067)	15000158	2
Unisan	Selenium (H1067)	15000158/0	4
Unisan	Selenium (H1067)	15000157	3
Unisan	Selenium (H1067)	65501/2	1
Unisan	Selenium (H1067)	65601/1	1
Unisan	Selenium (H1067)	65601	4
Unisan	Selenium (H1067)	65601/2	1
Unisan	Selenium (H1067)	65601	2
Unisan	Selenium (H1067)	15000159	1
Unisan	Selenium yeast (H1068)	16001450/0	1
Unisan	Selenium yeast (H1068)	25709	1
Unisan	Selenium yeast (H1068)	9001	1
Unisan	Silymarin (H1071)	950401	1
Unisan	Silymarin (H1071)	15001046	3
Unisan	Silymarin (H1071)	46101	1
Unisan	Silymarin (H1071)	15001046/1	1
Unisan	Silymarin (H1071)	67701/7	1
Unisan	Sodium citrate	S/WReferenzUmstell	1
Unisan	Sodium sulfate decahy...	11294605	1
Unisan	Taurine (H1084)	1294960201	2
Unisan	Taurine (H1084)	16000946/1	1
Unisan	Taurine (H1084)	28001	11
Unisan	Taurine (H1084)	14001767/0	1
Unisan	Threonine (L-) (H1052)	14001718/1	1
Unisan	Threonine (L-) (H1052)	14001718	2
Unisan	Threonine (L-) (H1052)	14001718/7	1
Unisan	Trace elements JK (H1...	15000525/0	1
Unisan	Trace elements JK (H1...	15000526	1
Unisan	Trace elements JK (H1...	1028956901	1
Unisan	Trace elements JK (H1...	11701	14
Unisan	Trace elements JK (H1...	1026956901	1
Unisan	Trace elements JK (H1...	17000870	2
Unisan	Trace elements JK (H1...	72901	1
Unisan	Trace elements SE (H1...	15000555	1
Unisan	Trace elements SE (H1...	15000553	1
Unisan	Trace elements SE (H1...	14000721/0	1
Unisan	Trace elements SE (H1...	14000722	2
Unisan	Trace elements SE (H1...	31301	7
Unisan	Trace elements SE (H1...	1024012401	1
Unisan	Trace elements SE (H1...	14000721	1
Unisan	Trace elements SE (H1...	15000556/0	1
Unisan	Trace elements SE (H1...	1024929202	1
Unisan	Tryptophan (L-) (H1053)	15000443	2
Unisan	Tryptophan (L-) (H1053)	950001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	Tryptophan (L-) (H1053)	74501	1
Unisan	Tryptophan (L-) (H1053)	15000443/0	2
Unisan	Tryptophan (L-) (H1053)	150000441	1
Unisan	Tryptophan (L-) (H1053)	74501/4	1
Unisan	Tryptophan (L-) (H1053)	27902	3
Unisan	Tryptophan (L-) (H1053)	74501	1
Unisan	Tryptophan (L-) (H1053)	1824950001	2
Unisan	Tyrosine (L-) (H1054)	74601/3	1
Unisan	Tyrosine (L-) (H1054)	29501	3
Unisan	Tyrosine (L-) (H1054)	979102	1
Unisan	Tyrosine (L-) (H1054)	29502	5
Unisan	Valine (L-) (H1056)	15001436/1	1
Unisan	Valine (L-) (H1056)	15001436	3
Unisan	Valine (L-) (H1056)	1807961101	2
Unisan	Valine (L-) (H1056)	25012	2
Unisan	Vitamin B1 (H1101)	16000698	1
Unisan	Vitamin B1 (H1101)	15001622	1
Unisan	Vitamin B1 (H1101)	15001622/0	3
Unisan	Vitamin B1 (H1101)	15001620/0	1
Unisan	Vitamin B1 (H1101)	14001112/0	1
Unisan	Vitamin B1 (H1101)	30202	2
Unisan	Vitamin B1 (H1101)	1622977301	1
Unisan	Vitamin B1 (H1101)	17001075	1
Unisan	Vitamin B1 (H1101)	977301	2
Unisan	Vitamin B12 (H1102)	75301/2	1
Unisan	Vitamin B12 (H1102)	75301	5
Unisan	Vitamin B12 (H1102)	3441377	1
Unisan	Vitamin B12 (H1102)	25017	9
Unisan	Vitamin B12 (H1102)	26401	1
Unisan	Vitamin B12 (H1102)	15000993/3	1
Unisan	Vitamin B12 (H1102)	15000993	4
Unisan	Vitamin B12 (H1102)	17000276	3
Unisan	Vitamin B12 (H1102)	1061025017	1
Unisan	Vitamin B6 (H1103)	20101	1
Unisan	Vitamin B6 (H1103)	14001576/2	2
Unisan	Vitamin B6 (H1103)	14001576	3
Unisan	Vitamin B6 (H1103)	15000279/2	1
Unisan	Vitamin C (H1104)	16000200	1
Unisan	Vitamin C (H1104)	16000198	2
Unisan	Vitamin C (H1104)	1085025020	4
Unisan	Vitamin C (H1104)	16000203	1
Unisan	Vitamin C (H1104)	25020	2
Unisan	Vitamin C (H1104)	44502	5
Unisan	Vitamin C (H1104)	16000231	1
Unisan	Vitamin C (H1104)	73901	6
Unisan	Vitamin C (H1104)	74001	2
Unisan	Vitamin C (H1104)	14001780/0	1
Unisan	Vitamin C (H1104)	14001788	1
Unisan	Vitamin C (H1104)	14001797	1
Unisan	Vitamin complex (H1100)	49302	3
Unisan	Vitamin complex (H1100)	14001590/1	1
Unisan	Vitamin complex (H1100)	1014009901	1
Unisan	Vitamin complex (H1100)	25201	1
Unisan	Vitamin complex (H1100)	25101	1
Unisan	Vitamin complex (H1100)	H110031	2
Unisan	Vitamin complex (H1100)	15001275/0	2
Unisan	Vitamin complex (H1100)	74101	1
Unisan	Vitamin complex (H1100)	16000962	2
Unisan	Vitamin complex (H1100)	17001750	2
Unisan	Vitamin complex (H1100)	77001	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan	Vitamin complex (H1100)	904006	1
Unisan	Vitamin complex (H1100)	77001	1
Unisan	Vitamin complex (H1100)	14001589	3
Unisan	Vitamin complex (H1100)	14001590/0	1
Unisan	Vitamin complex (H1100)	16000958/0	2
Unisan	Vitamin complex (H1100)	1400159/0	1
Unisan	Vitamin complex (H1100)	15001276/0	1
Unisan	Vitamin complex (H1100)	15001277/0	1
Unisan	Vitamin D3 (H1105)	16001212	3
Unisan	Vitamin D3 (H1105)	15000490/0	2
Unisan	Vitamin D3 (H1105)	14000481/2	3
Unisan	Vitamin D3 (H1105)	16001212/1	1
Unisan	Vitamin D3 (H1105)	16000952	2
Unisan	Vitamin D3 (H1105)	14000481	4
Unisan	Vitamin D3 (H1105)	1335008101	3
Unisan	Vitamin D3 (H1105)	1335008108	1
Unisan	Vitamin D3 (H1105)	1335959601	1
Unisan	Vitamin D3 (H1105)	8101	9
Unisan	Vitamin D3 (H1105)	15000490	4
Unisan	Vitamin E NAT (H1107)	16001246	2
Unisan	Vitamin E NAT (H1107)	16000591	2
Unisan	Vitamin E NAT (H1107)	71501	4
Unisan	Vitamin E NAT (H1107)	16000589/0	2
Unisan	Vitamin E NAT (H1107)	1094025021	1
Unisan	Vitamin E NAT (H1107)	1095961301	1
Unisan	Vitamin E NAT (H1107)	14001821	1
Unisan	Vitamin E NAT (H1107)	15001558	1
Unisan	Vitamin E NAT (H1107)	15001557/0	3
Unisan	Vitamin E NAT (H1107)	15001659/0	1
Unisan	Vitamin E NAT (H1107)	1094017701	2
Unisan	Vitamin K2-Ca-Mn comp...	69001	2
Unisan	Zinc (H1110)	1164991001	6
Unisan	Zinc (H1110)	16000196	5
Unisan	Zinc (H1110)	16000196/0	2
Unisan	Zinc (H1110)	991001	11
Unisan	Zinc (H1110)	14001753	1
Unisan	Zinc (H1110)	14001753/0	1
Unisan	Zinc (H1110)	14001803	3
Unisan	Zinc (H1110)	14001803/1	1
UNISAN	Antioxidants NAT comp...	14000042	1
UNISAN	Carnitine (L-) (H1038)	130605302	3
UNISAN	Carnitine (L-) (H1038)	130605300	6
UNISAN	Carnosine (L-) (H1014)	74801	1
UNISAN	Chromium (H1012)	14000429	6
UNISAN	Coenzyme Q10 (H1013)	14001546	4
UNISAN	Glucosamine sulfate (...)	34001	1
UNISAN	Glutathione (H1020)	74901	4
UNISAN	Iron (H1015)	14000285	2
UNISAN	Lipoic acid (a-) (H10...	14001726	2
UNISAN	Selenium (H1067)	65501	2
UNISAN	Vitamin B1 (H1101)	75001	1
UNISAN	Vitamin D3 (H1105)	63601	1
UNISAN	Vitamin E NAT (H1107)	71701	1
Unisan	Carnitine (L-) (H1038)	1184965901	1
Unisan	Copper (H1032)	1542933901	1
Unisan	Lysine (L-) (H1047)	1801028801	1
Unisan	Trace elements SE (H1...	1024042301	1
Unisan	Vitamin complex (H1100)	1014025401	1
Unisan	Zinc (H1110)	1164949701	1
Unisan 02.04.2016	Antioxidants NAT comp...	15000874/0	2

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Unisan 02.04.2016	Bioflavonoid complex ...	14001675/2	1
Unisan 02.04.2016	Coenzyme Q10 (H1013)	14001546/3	1
Unisan 02.04.2016	Isoflavones (H1025)	15000205/1	1
Unisan 02.04.2016	Vitamin E NAT (H1107)	15001558/0	1
Unisan 11.08.2105	Iron (H1015)	14000285/2	1
Unisan 2.10.2015	Bioflavonoid complex ...	14001675/1	2
Unisan 2.10.2015	Carnitine (L-) (H1038)	15000052/0	2
Unisan 2.10.2015	Coenzyme Q10 (H1013)	14001546/2	1
Unisan 2.10.2015	Iron (H1015)	14000285/3	2
Unisan 2.10.2015	Vitamin E NAT (H1107)	14001821/0	2
Unisan Gmbh, 78465 Ko...	Antioxidants NAT comp...	67301	3
Unisan Gmbh, 78465 Ko...	Arginine (L-) (H1036)	52401	4
Unisan Gmbh, 78465 Ko...	Creatine (H1031)	1830025003	1
Unisan Gmbh, 78465 Ko...	Folic acid (H1016)	1611030201	1
Unisan Gmbh, 78465 Ko...	Glutamine (L-) (H1043)	26904	2
Unisan Gmbh, 78465 Ko...	Vitamin D3 (H1105)	63701	9
Unisan/Hepart	Iron (H1015)	15001259/0	1
Unisan/Hepart	Vitamin D3 (H1105)	15000490/5	2
Unisan/Hepart AG	Aminomix NAC (H1002)	41402	7
Unisan/Hepart AG	Catechin extract (H10...	18601	4
Unisan/Hepart AG	Chaga (H1011)	46701	1
Unisan/Hepart AG	Chondroitin sulfate ...	31801	9
Unisan/Hepart AG	Ginseng (H1018)	38802	4
Unisan/Hepart AG	Glutamine (L-) (H1043)	26901	1
Unisan/Hepart AG	Glycine (H1021)	28101	7
Unisan/Hepart AG	Inositol hexanicotina...	46001	2
Unisan/Hepart AG	Isoflavones (H1025)	38801	2
Unisan/Hepart AG	Leucine (L-) (H1046)	52601	4
Unisan/Hepart AG	Lipoic acid (a-) (H10...	54301	2
Unisan/Hepart AG	Methylsulfonylmethane...	1819025013	1
Unisan/Hepart AG	Ornithine (L-) (H1048)	27001	2
Unisan/Hepart AG	Phenylalanine (L-) (H...	1828025009	1
Unisan/Hepart AG	Potassium (H1030)	39302	3
Unisan/Hepart AG	Silymarin (H1071)	29901	6
Unisan/Hepart AG	Vitamin B1 (H1101)	53301	1
Unisan/Hepart AG	Vitamin C (H1104)	44602	3
Unisan/Hepart AG	Vitamin complex (H1100)	49201	5
Unisan/Hepart AG	Vitamin complex (H1100)	48301	1
v	Progesterone, microni...	14253906	1
VDL	Basil, rubbed	12007004	1
VDL	Boldo leaf, cut	11016503	1
VDL	Chloramphenicol	140743-01	1
VDL	Citric acid monohydrate	263669	1
VDL	Clioquinol	131144131	1
VDL	Eucalyptus leaf, cut	2011B110801	4
VDL	Java tea, cut	2221B101213	1
VDL	Macrogols 4000	16308906	1
VDL	Primrose flower	11206806	1
VDL	Raspberry leaf, cut	10013406	3
VDL	Salicylic acid	11197122	1
VDL	Yellow chaste weed, cut	10244406	1
Vdl 09.05.14;Preis 7,45	Riboflavin	13305203/15R05	1
Vdl 14.07.14	Clotrimazole	13224225	3
Vdl 16.10.14 12,70 euro	Estradiol benzoate	14B20-02-292001-	1
Vdl 23.10.14 19,50 euro	Budesonide, micronized	14071707	3
Vdl 28.04.14;Preis 39...	Progesterone, microni...	13155902, 15R05	1
vdl 28.04.2014 preis:...	Sorbic acid	12330618, 15R05	1
Vdl 28.10.14 10,32 euro	Isosorbide dinitrate ...	14D09-B03-294280--	1
VDL;03.03.15;3,04EUR	Rice starch	14011627	1
VDL;15.06.15;6,83EUR	Caraway fruit	13306007	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
VDL;16.04.15;5,60EUR	Nicotinamide	1408019-01	3
VDL;23.6.14; 2.49EUR	Citric acid	13243213	21
VDL;5,78EUR;18.12	Sorbic acid	14176206	10
VEB Jenapharmn	Procaine hydrochloride	1011989	1
vedono	Aniseed, powdered	11306208	1
vedono	Frankincense	12120012	1
vedono	Frankincense	12120009	5
vedono	Hawthorn leaf with fl...	12225103	4
vedono	Rhubarb, cut	12211907	6
vedono	Roselle	12007313	1
vedono	Wild thyme, cut	11374807	2
vedono	Zedoary root, cut	13010106	3
Wachs neu	Wax, yellow		2
WE:05.04.2016 Phönix	Betamethasone dipropi...	15163304	12
WE:Sanacorp 02.09.16 ...	Mannitol	16212101	13
Weltecke	Angelica root	38011111	1
Weltecke	Meadowsweet herb, cut	67010413	1
Weltecke	Thyme, rubbed	67101212	1
Weltecke / GEHE	Primrose flower	63011112	1
Wep apot	Ascorbic acid		1
Wepa	Ascorbic acid	180316-C01	1
Wepa	Ascorbic acid	261120415	1
wepa Apoday	Ascorbic acid	2616110375431	1
Wepa Apothekenbedarf	Dequalinium chloride	151203202	1
Wepa Paulus	Ascorbic acid	11406201	1
wepa/Sanaco	Clioquinol	1439118	1
Werner	Sucrose	1901208731	2
Wolf	Alchemilla, cut	23187	4
Wolf	Alchemilla, cut	p931	1
Wolf	Alchemilla, cut	E-822	1
Wolf	Alchemilla, cut	B-698	2
Wolf	Aniseed	B-671	2
Wolf	Aniseed	L401F1413	1
Wolf	Aniseed	L401F292b	4
Wolf	Aniseed	B671	1
Wolf	Aniseed	215346	1
Wolf	Aniseed	B216346	14
Wolf	Aniseed	216346	1
Wolf	Bitter fennel	B492	1
Wolf	Bitter fennel	B-492	5
Wolf	Bitter fennel	B-362	1
Wolf	Bitter fennel	B-682	1
Wolf	Bitter fennel	B362	2
Wolf	Bitter fennel	B-493	1
Wolf	Bitter fennel	B382	1
Wolf	Black cumin seed		1
Wolf	Black cumin seed	207204	1
Wolf	Black cumin seed	20720	8
Wolf	Black cumin seed	B-471	1
Wolf	Black cumin seed	L402F1713	1
Wolf	Black cumin seed	L310F141b	6
Wolf	Black cumin seed	L402F171b	3
Wolf	Calendula flower	23342	1
Wolf	Caraway fruit	B-224	1
Wolf	Dill fruit	23622	6
Wolf	Dill fruit	L30057	5
Wolf	Dill fruit	E-501	1
Wolf	Elder flower, rubbed	B-810	1
Wolf	Elder flower, rubbed	22749	1
Wolf	Elder flower, rubbed	E-738	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
Wolf	Elder flower, rubbed	22765	1
Wolf	Equisetum stem, cut	B-695	1
Wolf	Equisetum stem, cut	B233885	4
Wolf	Lemon grass, cut	B654	1
Wolf	Lemon verbena leaf, cut	23635	1
Wolf	Lemon verbena leaf, cut	E-587	1
Wolf	Lemon verbena leaf, cut	B-312	1
Wolf	Lime flower, cut	B449	1
Wolf	Lime flower, cut	22094	1
Wolf	Marjoram herb, rubbed	P232243	9
Wolf	Marjoram herb, rubbed	163516	1
Wolf	Marjoram herb, rubbed	B163516	6
Wolf	Matricaria flower	23242	1
Wolf	Nettle leaf, cut	P234605	10
Wolf	Nettle leaf, cut	B821	1
Wolf	Nettle leaf, cut	B-821	1
Wolf	Peppermint leaf, cut	B-120	1
Wolf	Peppermint leaf, cut	B-725	2
Wolf	Peppermint leaf, cut	B725	1
Wolf	Peppermint leaf, cut	P221379	7
Wolf	Peppermint leaf, cut	B243	1
Wolf	Raspberry leaf, cut	P201727	19
Wolf	Raspberry leaf, cut	P222676	1
Wolf	Raspberry leaf, cut	B-457	1
Wolf	Sage leaf, cut	B735	1
Wolf	Shepherd's purse, cut	B449	2
Wolf	St. John's wort, cut	19468	2
Wolf	St. John's wort, cut	206152	1
Wolf	St. John's wort, cut	B-455	2
Wolf	St. John's wort, cut	B206152	1
Wolf	St. John's wort, cut	B-424	1
Wolf	Yarrow herb, cut	B230010	9
Wolf	Yarrow herb, cut	230010	1
Wolf Apotheke	Citric acid	58(18097609)	1
Wolf-Apotheke	Prednisolone acetate	14380911	1
x	Oxytetracycline hydro...		1
xxx	Camphor	709E-01221	1
xxx	Hydrocortisone	911Q-01719	2
xxx	Salicylic acid	4311A-01045	1
xxxxx	Mannitol	xxxx	1
Zentiva/Alliance Heal...	Ascorbic acid	15300007	8
-	Angelica root	130924	1
-	Bitter orange flower ...	2011A-01445	1
-	Burdock root, cut	5110I-02131	1
-	Chelidonium majus her...	810I-01632	1
-	Chickweed, cut	2011A-01686	1
-	Chicory herb, cut	4411I-01634	1
-	Elecampane rhizome, cut	712M-02181	1
-	Gentian root, powdered	130063	1
21.03.2017	Glucose	6910499	1
5N Plus Lübeck GmbH/ ...	Heavy alkaline bismut...	12146908	1
14026902	Hydroquinone	InterneCharge1214Q	1
-	Iceland moss, cut	4711E-01774	1
8i	Lactose monohydrate	309508/0	1
8i	Lactose monohydrate	309351/0	1
-	Meadowsweet herb, cut	1511M-01685	1
-	Meadowsweet herb, cut	130849	1
-	Melilot herb, cut	811Q-01664	1
-	Melissa leaf, cut	131146	1
?	Potassium sorbate	10012203	1

continued on the next page

continued from previous page

Supplier	Substance	Batch	Spectra
-	Raspberry leaf, cut	130849	1
14076202	Salicylic acid	704.2014	1
-	Silverweed, cut	711M-01619	1
-	Silverweed, cut	130849	1
-	Wild thyme, cut	130685	1
-	Wild thyme, cut	131023	1

Appendix D: Requirements of validation

In order to ensure adherence to the safe scientific status, the individual methods for manufacturing and testing must be validated under certain circumstances (compare § 34 para. 1 no. 3, § 35 para. 1 no. 4 and para. 4 sentence 1 no. 2 b, para. 6 sentence 3 *ApoBetrO* [Pharmacies Rules and Regulations]). The *ApoBetrO* [Pharmacies Rules and Regulations] incorporates a legal definition in § 1 a para. 16 (quotation translated):

“Validation is the provision of documented proof which with a high degree of safety documents that, via a specific process or standard work process, a medicinal product is manufactured and tested, which is in accordance with previously determined quality features.”

Validation documentation can be used to prove that methods or devices which are not described in the Pharmacopoeia within the meaning of § 6 para. 1 sentence 3 *ApBetrO* [Pharmacies Rules and Regulations] achieve the same results as those in the Pharmacopoeia. On the other hand, with the requirements of the demanded validation it must be observed whether the respective testing method is already incorporated in the Pharmacopoeia.

NIR spectroscopy as a general testing method need not be validated in accordance with the express ruling in the *Ph. Eur. Section 1.1* [3], as it is already described in *Section 2.2.40* of the *Ph. Eur.* as an area of application for the identification of raw materials.

However, a special validation requirement exists for the reference database. This requirement is met with the existing document. Further requirements or rules as to how this proof must be furnished do not exist. It is required that the processes guarantee the same results as the methods and devices in the Pharmacopoeia [17].

Carrying out identity tests with *Apo-Ident* is therefore also possible if the NIR spectroscopy process is not required in the Pharmacopoeia monograph of the substance for identity testing. All NIR analyses with *Apo-Ident* prove several, often all molecule groups and are therefore comparable with a series of individual, targeted chemical proofs [4]. Therefore, the identity proof with *Apo-Ident* replaces the monograph test series (with two or more test combinations).

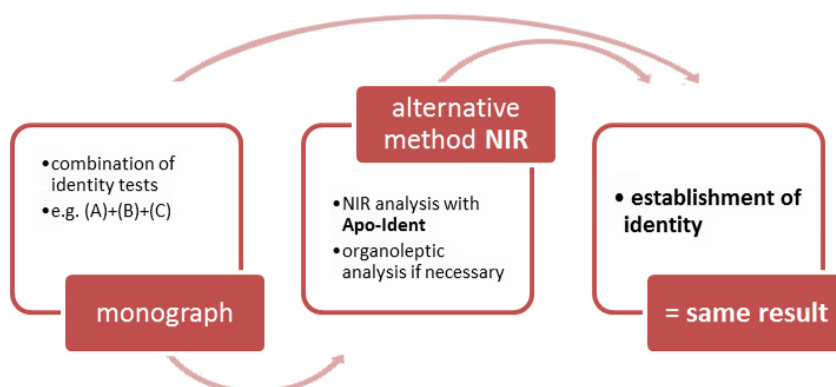


Figure 2: The combination of tests of the monograph is replaced by the alternative method NIR spectroscopy using *Apo-Ident*. This is permissible because both test procedures result in the establishment of the identity of the raw material.

With this validation documentation, proof is furnished that identical results are achieved with *Apo-Ident* and Pharmacopoeia methods, i.e. confirmation of the identity of the raw material [2].

Appendix E: Conformity of Apo-Ident with the European pharmacopeia

According to *Ph. Eur. Section 2.2.40*, NIR spectroscopy is basically suitable for: “Identification of agents, excipients, dosage forms, intermediate manufacturing products, chemical raw materials and packaging materials” ([3], quotation translated).

The fact that *Apo-Ident* meets the further criteria of the European Pharmacopoeia under the headings in *Section 2.2.40*

- Apparatus
- Measurement methods
- Sample preparation and presentation
- Testing the functionality of the instrument
- Identification and characterisation (qualitative analysis)
- Quantitative analysis
- Ongoing model evaluation
- Transfer of databases
- Data storage

can be proven based on the *HiperScan GmbH* documentation of “Meeting 2.2.40 *Ph. Eur.* by *Apo-Ident*” [4].

Appendix F: Addition to the second-stage models

As described in the general part of the *Apo-Ident* validation documentation, in the section [Model creation procedure and validation runs](#), under 2. *Generating the chemometric models (calibration) in f)*, it is possible to create a further chemometric model for certain substances (second-stage model) and to undertake the assessment in several steps.

Substances in second-stage models

In the *Apo-Ident Update 2019-05*, a distinction is made between the following substances in the listed second-stage models:

Submodel 1

Bacitracin
Betamethasone, micronized
Betamethasone valerate
Budesonide, micronized
Capsaicin
Dexamethasone
Erythromycin
Heparin sodium
Norethisterone acetate
Prednicarbate, micronized
Prednisolone, micronized
Prednisone
Triamcinolone, micronized
Triamcinolone acetonide

Submodel 2

Beclometasone dipropionate, anhydrous
Betamethasone dipropionate, micronized
Clobetasol propionate
Diphenylcyclopropenone
Gentamicin sulfate
Sodium benzoate
Sodium citrate
Oxybutynin hydrochloride

Submodel 3

Alfatradiol
Estriol
Pregnenolone
Progesterone, micronized
Quinine hydrochloride
Quinine sulfate dihydrate
Spironolactone
Testosterone propionate

Liquid submodel

2-ethylhexyl laurate
Base cream DAC

Eucerin® cum Aqua
Lanette® ointment (conserved)
Non-ionic hydrophil ointment
Squalane
Wool alcohols cream DAB
Wool alcohols cream SR

References

- [1] ABDA – BUNDESVEREINIGUNG DEUTSCHER APOTHEKERVERBÄNDE: Verordnung über den Betrieb von Apotheken (Apothekenbetriebsordnung – ApBetrO), 2012
- [2] REIMANN, B. ; REGIERUNGSPRÄSIDIUM DARMSTADT: Hinweise zur ordnungsgemäßen Prüfung von Arzneimitteln und Ausgangsstoffen (§§ 6 und 11 *ApBetrO*), 2007
- [3] *Europäisches Arzneibuch, Grundwerk 2017*. 9. Ausgabe. Deutscher Apotheker Verlag (978-3-7692-6679-5)
- [4] HIPERSCAN GMBH: Erfüllung von 2.2.40 Ph. Eur. durch Apo-Ident, 2013
- [5] ARBEITSGEMEINSCHAFT DER PHARMAZIERÄTE DEUTSCHLANDS (APD): Resolution in Verbindung mit §§ 6 und 11 *ApBetrO* Verwendung eines Nah-Infrarot-Spektrometers (NIR) zur Identitätsprüfung, 16. 10. 2013, DAZ 48, November 2013
- [6] ARBEITSGEMEINSCHAFT DER PHARMAZIERÄTE DEUTSCHLANDS (APD): Resolution 2014, Arbeitsgemeinschaft der Pharmazieräte Deutschlands (APD), Oktober 2014
- [7] *DAC/NRF*. Govi-Verlag (978-3-7741-0044-2)
- [8] KESSLER, W.: *Multivariate Datenanalyse*. WILEY-VCH Verlag, 2007 (978-3-527-31262-7)
- [9] NÆS, T. ; ISAKSSON, T. ; FEARN, T. ; DAVIES, T.: *Multivariate Calibration and Classification*. NIR Publications, 2002 (978 0 9528666 2 6)
- [10] HANLEY, J. A. ; LIPPMAN-HAND, A.: If nothing goes wrong, is everything all right? In: *Journal of the American Medical Association* 249 (1983), S. 1743–1745
- [11] JOVANOVIĆ, B. D. ; LEVY, P. S.: A Look at the Rule of Three. In: *The American Statistician* 51 (1997), S. 137–139
- [12] BRONSTEIN, I. N. ; SEMENDJAJEW, K. A. ; MUSIOL, G. ; MÜHLIG, H.: *Taschenbuch der Mathematik*. 5. überarbeitete und erweiterte Auflage. Verlag Harri Deutsch, 2000 (3-8171-2015-2)
- [13] MAHALANOBIS, P.: On the generalized distance in statistics. In: *Proc. Nat. Inst. Sci. India (Calcutta)* 2 (1936), S. 49–55
- [14] YAMBOR, B. ; DRAPER, W. ; BEVERIDGE, R.: Analyzing PCA-based face recognition algorithms: Eigenvector selection and distance measures. In: *Second Workshop Empirical Evaluation in Computer Vision* (2000)
- [15] HIPERSCAN GMBH: Identifikationsmethodik Apo-Ident, 2012
- [16] HIPERSCAN GMBH: Datenvorbehandlung des Identifikationssystems Apo-Ident, 2012
- [17] CYRAN ; RÖTTA: Apothekenbetriebsordnung, Kommentar § 6 Rn. 10, 2010

Index

- Acidum phenylaethylbarbituricum, *see* [Phenobarbital](#)
- Amfetamini sulfas, *see* [Amphetamine sulfate](#)
- Amphetamine sulfate, [19](#)
- Amphetamine sulphate, *see* [Amphetamine sulfate](#)
- Barbital sodium, [25](#)
- Barbitalum natriucum, *see* [Barbital sodium](#)
- Cocaine hydrochloride, [29](#)
- Cocainum hydrochloricum, *see* [Cocaine hydrochloride](#)
- Codeine phosphate hemihydrate, [33](#)
- Codeini phosphas hemihydricus, *see* [Codeine phosphate hemihydrate](#)
- Dexamphetamine sulfate, [37](#)
- Dexamphetamine sulphate, *see* [Dexamphetamine sulfate](#)
- Dexamphetaminum sulfuricum, *see* [Dexamphetamine sulfate](#)
- Dihydrocodeine hydrogen tartrate, [41](#)
- Dihydrocodeini hydrogenotartras, *see* [Dihydrocodeine hydrogen tartrate](#)
- Fentanili citras, *see* [Fentanyl citrate](#)
- Fentanyl citrate, [45](#)
- Hydromorphinonum hydrochloricum, *see* [Hydromorphone hydrochloride](#)
- Hydromorphone hydrochloride, [49](#)
- Methadone hydrochloride, [53](#)
- Methadoni hydrochloridum, *see* [Methadone hydrochloride](#)
- Midazolam, [61](#)
- Midazolam hydrochloride, [65](#)
- Midazolami hydrochloridum, *see* [Midazolam hydrochloride](#)
- Midazolamum, *see* [Midazolam](#)
- Morphine hydrochloride, [69](#)
- Morphine sulfate, [75](#)
- Morphine sulphate, *see* [Morphine sulfate](#)
- Morphini hydrochloridum, *see* [Morphine hydrochloride](#)
- Morphini sulfas, *see* [Morphine sulfate](#)
- Oxycodone hydrochloride, [79](#)
- Oxycodoni hydrochloridum, *see* [Oxycodone hydrochloride](#)
- Phenobarbital, [83](#)